

# GDCM

## 2.4.4

Generated by Doxygen 1.8.9.1

Sat Sep 26 2015 07:03:31



# Contents

<b>1</b>	<b>GDCM Documentation</b>	<b>1</b>
<b>2</b>	<b>off-screen rendering of DICOM images</b>	<b>3</b>
2.1	SYNOPSIS . . . . .	3
2.2	DESCRIPTION . . . . .	3
2.3	PARAMETERS . . . . .	3
2.4	OPTIONS . . . . .	3
2.4.1	OPTIONS . . . . .	3
2.4.2	general options . . . . .	3
2.5	Simple usage . . . . .	4
2.6	SEE ALSO . . . . .	4
2.7	COPYRIGHT . . . . .	4
<b>3</b>	<b>Convert a file supported by VTK into DICOM.</b>	<b>5</b>
3.1	SYNOPSIS . . . . .	5
3.2	DESCRIPTION . . . . .	5
3.3	PARAMETERS . . . . .	5
3.4	OPTIONS . . . . .	5
3.4.1	OPTIONS . . . . .	5
3.4.2	compression options . . . . .	6
3.4.3	general options . . . . .	6
3.4.4	environment variable . . . . .	6
3.5	DESCRIPTION . . . . .	6
3.5.1	CONVERT Metalmage (mhd, mha) . . . . .	6
3.5.2	CONVERT MHA/MHD . . . . .	7
3.5.3	CONVERT VTI . . . . .	7
3.5.4	CONVERT VTK . . . . .	7
3.6	CONVERT DICOM . . . . .	7
3.7	RoundTrip DICOM to MHD to DICOM . . . . .	7

3.8	gdcm2vtk notes	7
3.9	SEE ALSO	8
3.10	COPYRIGHT	8
<b>4</b>	<b>Tool to anonymize a DICOM file.</b>	<b>9</b>
4.1	SYNOPSIS	9
4.2	DESCRIPTION	9
4.3	PARAMETERS	9
4.4	OPTIONS	10
4.4.1	Required parameters	10
4.4.2	OPTIONS	10
4.4.3	encryption options	10
4.4.4	dumb mode options	10
4.4.5	general options	10
4.4.6	environment variable	11
4.5	Typical usage	11
4.5.1	De-identification (anonymization, encrypt)	11
4.5.2	Re-identification (de-anonymization,decrypt)	11
4.5.3	Multiple files caveat	11
4.5.4	Dumb mode	11
4.5.4.1	Irreversible Anonymization	12
4.6	OpenSSL	12
4.6.1	Generating a Private Key	12
4.6.2	Generating a Certificate	13
4.7	DICOM Standard:	13
4.8	Warnings	13
4.9	SEE ALSO	13
4.10	COPYRIGHT	13
<b>5</b>	<b>Tool to convert DICOM to DICOM.</b>	<b>15</b>
5.1	SYNOPSIS	15
5.2	DESCRIPTION	15
5.3	PARAMETERS	15
5.4	OPTIONS	15
5.4.1	PARAMETERS	15
5.4.2	OPTIONS	15
5.4.3	image options	16
5.4.4	JPEG options	16



5.4.5	JPEG-LS options	16
5.4.6	J2K options	16
5.4.7	general options	16
5.4.8	special options	16
5.4.9	environment variable	17
5.5	Simple usage	17
5.6	Typical usage	17
5.6.1	File Meta Header	17
5.6.2	Conversion to Explicit Transfer Syntax	18
5.6.3	Compressing to lossless JPEG	18
5.6.4	Compressing to lossy JPEG	18
5.6.5	Compressing to lossless JPEG-LS	18
5.6.6	Compressing to lossy JPEG-LS	18
5.6.7	Compressing to lossless J2K	18
5.6.8	Compressing to lossy J2K	18
5.6.9	Compressing to lossless RLE	19
5.6.10	Split encapsulated DICOM:	19
5.6.11	Forcing (re)compression	19
5.6.12	Decompressing a Compressed DICOM	19
5.6.13	Compressing an uncompressed Icon	19
5.6.14	Generating an Icon	20
5.6.15	Changing the planar Configuration	20
5.7	Lossless Conversion	20
5.8	Quality Control	20
5.8.1	DCMTK / dicom3tools	20
5.8.2	VIM: vimdiff	21
5.8.3	vbindiff	21
5.9	SEE ALSO	21
5.10	COPYRIGHT	21
<b>6</b>	<b>dumps differences of two DICOM files</b>	<b>23</b>
6.1	SYNOPSIS	23
6.2	DESCRIPTION	23
6.3	PARAMETERS	23
6.4	OPTIONS	23
6.4.1	OPTIONS	23
6.4.2	general options	23

6.5	Simple usage	24
6.6	SEE ALSO	24
6.7	COPYRIGHT	24
<b>7</b>	<b>dumps a DICOM file, it will display the structure and values contained in the specified DICOM file.</b>	<b>25</b>
7.1	SYNOPSIS	25
7.2	DESCRIPTION	25
7.3	PARAMETERS	25
7.4	OPTIONS	25
7.4.1	OPTIONS	25
7.4.2	general options	26
7.4.3	special options	26
7.5	Typical usage	26
7.5.1	Printing Implicit Transfer Syntax	26
7.5.2	Print Private Attributes	27
7.5.3	SIEMENS CSA Header	27
7.5.4	GEMS Protocol Data Block	27
7.5.5	ELSCINT Protocol Information	28
7.5.6	VEPRO Protocol Information	28
7.5.7	Philips Private MR Series Data Storage (1.3.46.670589.11.0.0.12.2)	29
7.5.8	Encapsulated ASN1 Structure	30
7.6	SEE ALSO	31
7.7	COPYRIGHT	31
<b>8</b>	<b>Tool to generate a DICOMDIR file from a File-Set.</b>	<b>33</b>
8.1	SYNOPSIS	33
8.2	DESCRIPTION	33
8.3	PARAMETERS	33
8.4	OPTIONS	33
8.4.1	Parameters	33
8.4.2	OPTIONS	33
8.4.3	general options	33
8.4.4	environment variable	34
8.5	Typical usage	34
8.6	NOTE	34
8.7	SEE ALSO	34
8.8	COPYRIGHT	34

<b>9</b>	<b>Manipulate DICOM image file.</b>	<b>35</b>
9.1	SYNOPSIS	35
9.2	DESCRIPTION	35
9.3	PARAMETERS	35
9.4	OPTIONS	35
9.4.1	PARAMETERS	35
9.4.2	OPTIONS	36
9.4.3	fill options	36
9.4.4	general options	36
9.4.5	environment variable	36
9.5	Supported File Format (appropriate file extension) <code>gdcmimg</code>	36
9.6	Typical usage	37
9.6.1	Remove a rectangular part of the image	37
9.6.2	Convert RAW to DICOM	37
9.6.3	Convert PGM/PNM/PPM to DICOM	38
9.6.4	Convert RLE to DICOM	38
9.6.5	Convert JPEG to DICOM	38
9.6.6	Convert J2K to DICOM	38
9.6.7	Specifying a SOP Class UID	38
9.7	Multiple Files	38
9.8	Start Offset	39
9.9	Warning	39
9.10	SEE ALSO	39
9.11	COPYRIGHT	39
<b>10</b>	<b>Display meta info about the input DICOM file.</b>	<b>41</b>
10.1	SYNOPSIS	41
10.2	DESCRIPTION	41
10.3	PARAMETERS	41
10.4	OPTIONS	41
10.4.1	OPTIONS	41
10.4.2	general options	41
10.4.3	environment variable	42
10.5	Simple usage	42
10.5.1	<code>gdcmData</code>	42
10.5.2	Davie Clunie datasets:	42
10.5.3	Checking the md5sum of the Pixel Data	43

10.5.4 Checking if Pixel Data is lossless . . . . .	43
10.6 SEE ALSO . . . . .	43
10.7 COPYRIGHT . . . . .	43
<b>11 Tool to convert PYPYRUS 3.0 to DICOM.</b>	<b>45</b>
11.1 SYNOPSIS . . . . .	45
11.2 DESCRIPTION . . . . .	45
11.3 PARAMETERS . . . . .	45
11.4 OPTIONS . . . . .	45
11.4.1 PARAMETERS . . . . .	45
11.4.2 OPTIONS . . . . .	45
11.4.3 general options . . . . .	45
11.4.4 environment variable . . . . .	46
11.5 Simple usage . . . . .	46
11.6 SEE ALSO . . . . .	46
11.7 COPYRIGHT . . . . .	46
<b>12 Tool to convert PDF to PDF/DICOM.</b>	<b>47</b>
12.1 SYNOPSIS . . . . .	47
12.2 DESCRIPTION . . . . .	47
12.3 PARAMETERS . . . . .	47
12.4 OPTIONS . . . . .	47
12.4.1 general options . . . . .	47
12.5 Usage Example . . . . .	48
12.6 PDF Info Mapping . . . . .	48
12.7 SEE ALSO . . . . .	49
12.8 COPYRIGHT . . . . .	49
<b>13 Extract Data Element Value Field.</b>	<b>51</b>
13.1 SYNOPSIS . . . . .	51
13.2 DESCRIPTION . . . . .	51
13.3 PARAMETERS . . . . .	51
13.4 OPTIONS . . . . .	51
13.4.1 PARAMETERS . . . . .	51
13.4.2 OPTIONS . . . . .	51
13.4.3 general options . . . . .	51
13.5 Typical usage . . . . .	52
13.5.1 Copy Attribute Value to file . . . . .	52

13.5.2 Extract Pixel Data . . . . .	52
13.5.3 Encapsulated Syntax . . . . .	52
13.5.4 Extract fragments as single file . . . . .	53
13.6 Footnote about JPEG files . . . . .	54
13.7 SEE ALSO . . . . .	54
13.8 COPYRIGHT . . . . .	54
<b>14 Scan a directory containing DICOM files.</b>	<b>55</b>
14.1 SYNOPSIS . . . . .	55
14.2 DESCRIPTION . . . . .	55
14.2.1 PARAMETERS . . . . .	55
14.2.2 OPTIONS . . . . .	55
14.2.3 general options . . . . .	55
14.3 Typical usage . . . . .	56
14.4 Simple usage . . . . .	56
14.5 Complex usage . . . . .	56
14.6 SEE ALSO . . . . .	56
14.7 COPYRIGHT . . . . .	56
<b>15 Tool to execute a DICOM Query/Retrieve operation</b>	<b>57</b>
15.1 SYNOPSIS . . . . .	57
15.2 DESCRIPTION . . . . .	57
15.3 PARAMETERS . . . . .	57
15.4 OPTIONS . . . . .	57
15.4.1 OPTIONS . . . . .	57
15.4.2 mode options . . . . .	57
15.4.3 C-STORE options . . . . .	58
15.4.4 C-FIND/C-MOVE options . . . . .	58
15.4.5 C-MOVE options . . . . .	58
15.4.6 general options . . . . .	58
15.4.7 environment variable . . . . .	58
15.5 C-ECHO usage . . . . .	59
15.6 C-STORE usage . . . . .	59
15.7 C-FIND usage . . . . .	59
15.8 C-MOVE usage . . . . .	60
15.9 patientroot notes . . . . .	60
15.10 Debugging . . . . .	60
15.11 Port Warning . . . . .	60

15.12C-STORE Warnings . . . . .	61
15.13C-MOVE Warnings . . . . .	61
15.14C-FIND IMAGE level (Composite Object Instance) . . . . .	61
15.15Storing the Query . . . . .	61
15.16DICOM Public Servers . . . . .	62
15.17SEE ALSO . . . . .	62
15.18COPYRIGHT . . . . .	62
<b>16 Concatenate/Extract DICOM files. . . . .</b>	<b>63</b>
16.1 SYNOPSIS . . . . .	63
16.2 DESCRIPTION . . . . .	63
16.3 PARAMETERS . . . . .	63
16.4 OPTIONS . . . . .	63
16.4.1 OPTIONS . . . . .	63
16.4.2 general options . . . . .	63
16.4.3 environment variable . . . . .	64
16.5 Typical usage . . . . .	64
16.5.1 SIEMENS Mosaic . . . . .	64
16.6 SEE ALSO . . . . .	65
16.7 COPYRIGHT . . . . .	65
<b>17 Simple DICOM viewer. . . . .</b>	<b>67</b>
17.1 SYNOPSIS . . . . .	67
17.2 DESCRIPTION . . . . .	67
17.3 PARAMETERS . . . . .	67
17.4 OPTIONS . . . . .	67
17.4.1 OPTIONS . . . . .	67
17.4.2 general options . . . . .	67
17.5 Typical usage . . . . .	68
17.6 Simple usage . . . . .	68
17.7 Wiki Link . . . . .	68
17.8 SEE ALSO . . . . .	68
17.9 COPYRIGHT . . . . .	68
<b>18 provides a tool to convert a DICOM file into a XML infoset and vice-versa. . . . .</b>	<b>69</b>
18.1 SYNOPSIS . . . . .	69
18.2 DESCRIPTION . . . . .	69
18.3 PARAMETERS . . . . .	69

18.4	OPTIONS	69
18.4.1	PARAMETERS	69
18.4.2	Options for DICOM to XML:	69
18.4.3	Options for XML to DICOM:	70
18.4.4	general options	70
18.5	SEE ALSO	70
18.6	COPYRIGHT	70
<b>19</b>	<b>Todo List</b>	<b>71</b>
<b>20</b>	<b>Deprecated List</b>	<b>73</b>
<b>21</b>	<b>Bug List</b>	<b>75</b>
<b>22</b>	<b>Namespace Index</b>	<b>77</b>
22.1	Namespace List	77
<b>23</b>	<b>Hierarchical Index</b>	<b>79</b>
23.1	Class Hierarchy	79
<b>24</b>	<b>Class Index</b>	<b>87</b>
24.1	Class List	87
<b>25</b>	<b>File Index</b>	<b>103</b>
25.1	File List	103
<b>26</b>	<b>Namespace Documentation</b>	<b>109</b>
26.1	gdcm Namespace Reference	109
26.1.1	Detailed Description	123
26.1.2	Typedef Documentation	124
26.1.2.1	AEComp	124
26.1.2.2	ASComp	124
26.1.2.3	BOOL_FUNCTION_PFILE_PFILE_POINTER	124
26.1.2.4	CSComp	124
26.1.2.5	DAComp	124
26.1.2.6	DTComp	124
26.1.2.7	FileList	124
26.1.2.8	IconImage	124
26.1.2.9	LOComp	124
26.1.2.10	LTComp	124

26.1.2.11 MacroEntry . . . . .	124
26.1.2.12 NestedMacroEntries . . . . .	124
26.1.2.13 PNComp . . . . .	124
26.1.2.14 SHComp . . . . .	124
26.1.2.15 STComp . . . . .	124
26.1.2.16 TMComp . . . . .	124
26.1.2.17 UIComp . . . . .	124
26.1.2.18 UTComp . . . . .	124
26.1.3 Enumeration Type Documentation . . . . .	124
26.1.3.1 CompOperators . . . . .	124
26.1.3.2 ECharSet . . . . .	125
26.1.3.3 EQueryLevel . . . . .	125
26.1.3.4 EQueryType . . . . .	125
26.1.3.5 ERootType . . . . .	126
26.1.3.6 LodModeType . . . . .	126
26.1.4 Function Documentation . . . . .	126
26.1.4.1 backslash . . . . .	126
26.1.4.2 GetVRFromTag . . . . .	126
26.1.4.3 operator"!=" . . . . .	126
26.1.4.4 operator"!=" . . . . .	126
26.1.4.5 operator<< . . . . .	126
26.1.4.6 operator<< . . . . .	126
26.1.4.7 operator<< . . . . .	126
26.1.4.8 operator<< . . . . .	126
26.1.4.9 operator<< . . . . .	126
26.1.4.10 operator<< . . . . .	127
26.1.4.11 operator<< . . . . .	127
26.1.4.12 operator<< . . . . .	127
26.1.4.13 operator<< . . . . .	127
26.1.4.14 operator<< . . . . .	127
26.1.4.15 operator<< . . . . .	127
26.1.4.16 operator<< . . . . .	127
26.1.4.17 operator<< . . . . .	127
26.1.4.18 operator<< . . . . .	127
26.1.4.19 operator<< . . . . .	127
26.1.4.20 operator<< . . . . .	127
26.1.4.21 operator<< . . . . .	127



26.1.4.22 operator<<	127
26.1.4.23 operator<<	127
26.1.4.24 operator<<	128
26.1.4.25 operator<<	128
26.1.4.26 operator<<	128
26.1.4.27 operator<<	128
26.1.4.28 operator<<	128
26.1.4.29 operator<<	128
26.1.4.30 operator<<	128
26.1.4.31 operator<<	128
26.1.4.32 operator<<	128
26.1.4.33 operator<<	128
26.1.4.34 operator<<	128
26.1.4.35 operator<<	128
26.1.4.36 operator<<	128
26.1.4.37 operator<<	128
26.1.4.38 operator<<	128
26.1.4.39 operator<<	128
26.1.4.40 operator<<	128
26.1.4.41 operator<<	129
26.1.4.42 operator<<	129
26.1.4.43 operator<<	129
26.1.4.44 operator<<	129
26.1.4.45 operator<<	129
26.1.4.46 operator<<	129
26.1.4.47 operator<<	129
26.1.4.48 operator<<	129
26.1.4.49 operator<<	129
26.1.4.50 operator<<	129
26.1.4.51 operator<<	129
26.1.4.52 operator<<	130
26.1.4.53 operator<<	130
26.1.4.54 operator<<	130
26.1.4.55 operator<<	130
26.1.4.56 operator<<	130
26.1.4.57 operator<<	130
26.1.4.58 operator<<	130

26.1.4.59 operator==	130
26.1.4.60 operator>>	130
26.1.4.61 operator>>	130
26.1.4.62 operator>>	130
26.1.4.63 to_string	130
26.1.4.64 TYPETOENCODING	131
26.1.5 Variable Documentation	131
26.1.5.1 GlobalInstance	131
26.1.5.2 VRINARY	131
26.2 gdcm::network Namespace Reference	131
26.2.1 Enumeration Type Documentation	135
26.2.1.1 EEventID	135
26.2.1.2 EStateID	135
26.2.2 Function Documentation	136
26.2.2.1 GetStateIndex	136
26.2.3 Variable Documentation	136
26.2.3.1 cMaxEventID	136
26.2.3.2 cMaxStateID	136
26.3 gdcm::SegmentHelper Namespace Reference	136
26.4 gdcm::terminal Namespace Reference	136
26.4.1 Detailed Description	137
26.4.2 Enumeration Type Documentation	137
26.4.2.1 Attribute	137
26.4.2.2 Color	137
26.4.2.3 Mode	138
26.4.3 Function Documentation	138
26.4.3.1 setattribute	138
26.4.3.2 setbgcolor	138
26.4.3.3 setfgcolor	138
26.4.3.4 setmode	138
<b>27 Class Documentation</b>	<b>139</b>
27.1 gdcm::network::AAbortPDU Class Reference	139
27.1.1 Detailed Description	140
27.1.2 Constructor & Destructor Documentation	140
27.1.2.1 AAbortPDU	140
27.1.3 Member Function Documentation	140

27.1.3.1	IsLastFragment	140
27.1.3.2	Print	140
27.1.3.3	Read	140
27.1.3.4	SetReason	141
27.1.3.5	SetSource	141
27.1.3.6	Size	141
27.1.3.7	Write	141
27.2	gdcm::network::AAssociateACPDU Class Reference	141
27.2.1	Detailed Description	142
27.2.2	Member Typedef Documentation	143
27.2.2.1	SizeType	143
27.2.3	Constructor & Destructor Documentation	143
27.2.3.1	AAssociateACPDU	143
27.2.4	Member Function Documentation	143
27.2.4.1	AddPresentationContextAC	143
27.2.4.2	GetNumberOfPresentationContextAC	143
27.2.4.3	GetPresentationContextAC	143
27.2.4.4	GetUserInfo	143
27.2.4.5	InitFromRQ	143
27.2.4.6	IsLastFragment	143
27.2.4.7	Print	143
27.2.4.8	Read	143
27.2.4.9	SetCalledAETitle	143
27.2.4.10	SetCallingAETitle	143
27.2.4.11	Size	143
27.2.4.12	Write	143
27.2.5	Friends And Related Function Documentation	143
27.2.5.1	AAssociateRQPDU	144
27.3	gdcm::network::AAssociateRJPDU Class Reference	144
27.3.1	Detailed Description	145
27.3.2	Constructor & Destructor Documentation	145
27.3.2.1	AAssociateRJPDU	145
27.3.3	Member Function Documentation	145
27.3.3.1	IsLastFragment	145
27.3.3.2	Print	145
27.3.3.3	Read	145
27.3.3.4	Size	145

27.3.3.5	Write	145
27.4	gdcm::network::AAssociateRQPDU Class Reference	145
27.4.1	Detailed Description	147
27.4.2	Member Typedef Documentation	147
27.4.2.1	PresentationContextArrayType	147
27.4.2.2	SizeType	147
27.4.3	Constructor & Destructor Documentation	147
27.4.3.1	AAssociateRQPDU	147
27.4.3.2	AAssociateRQPDU	147
27.4.4	Member Function Documentation	147
27.4.4.1	AddPresentationContext	147
27.4.4.2	GetCalledAETitle	148
27.4.4.3	GetCallingAETitle	148
27.4.4.4	GetNumberOfPresentationContext	148
27.4.4.5	GetPresentationContext	148
27.4.4.6	GetPresentationContextByAbstractSyntax	148
27.4.4.7	GetPresentationContextByID	148
27.4.4.8	GetPresentationContexts	148
27.4.4.9	GetReserved43_74	148
27.4.4.10	GetUserInformation	148
27.4.4.11	IsAETitleValid	148
27.4.4.12	IsLastFragment	148
27.4.4.13	Print	148
27.4.4.14	Read	148
27.4.4.15	SetCalledAETitle	148
27.4.4.16	SetCallingAETitle	148
27.4.4.17	SetUserInformation	149
27.4.4.18	Size	149
27.4.4.19	Write	149
27.4.5	Friends And Related Function Documentation	149
27.4.5.1	AAssociateACPDU	149
27.5	gdcm::AbortEvent Class Reference	149
27.6	gdcm::network::AbstractSyntax Class Reference	150
27.6.1	Detailed Description	150
27.6.2	Constructor & Destructor Documentation	151
27.6.2.1	AbstractSyntax	151
27.6.3	Member Function Documentation	151

27.6.3.1	GetAsDataElement	151
27.6.3.2	GetName	151
27.6.3.3	operator==	151
27.6.3.4	Print	151
27.6.3.5	Read	151
27.6.3.6	SetName	151
27.6.3.7	SetNameFromUID	151
27.6.3.8	Size	151
27.6.3.9	Write	151
27.7	gdcm::AnonymizeEvent Class Reference	151
27.7.1	Detailed Description	153
27.7.2	Member Typedef Documentation	153
27.7.2.1	Self	153
27.7.2.2	Superclass	153
27.7.3	Constructor & Destructor Documentation	153
27.7.3.1	AnonymizeEvent	153
27.7.3.2	~AnonymizeEvent	153
27.7.3.3	AnonymizeEvent	153
27.7.4	Member Function Documentation	153
27.7.4.1	CheckEvent	153
27.7.4.2	GetEventName	153
27.7.4.3	GetTag	153
27.7.4.4	MakeObject	153
27.7.4.5	SetTag	153
27.8	gdcm::Anonymizer Class Reference	154
27.8.1	Detailed Description	155
27.8.2	Constructor & Destructor Documentation	156
27.8.2.1	Anonymizer	156
27.8.2.2	~Anonymizer	156
27.8.3	Member Function Documentation	156
27.8.3.1	BALCPPProtect	156
27.8.3.2	BasicApplicationLevelConfidentialityProfile	156
27.8.3.3	CanEmptyTag	157
27.8.3.4	ClearInternalUIDs	157
27.8.3.5	Empty	157
27.8.3.6	GetBasicApplicationLevelConfidentialityProfileAttributes	157
27.8.3.7	GetCryptographicMessageSyntax	157

27.8.3.8	GetFile	157
27.8.3.9	New	157
27.8.3.10	RecurseDataSet	157
27.8.3.11	Remove	157
27.8.3.12	RemoveGroupLength	157
27.8.3.13	RemovePrivateTags	158
27.8.3.14	RemoveRetired	158
27.8.3.15	Replace	158
27.8.3.16	Replace	158
27.8.3.17	SetCryptographicMessageSyntax	158
27.8.3.18	SetFile	158
27.9	gdcm::AnyEvent Class Reference	158
27.10	gdcm::network::ApplicationContext Class Reference	160
27.10.1	Detailed Description	160
27.10.2	Constructor & Destructor Documentation	160
27.10.2.1	ApplicationContext	160
27.10.3	Member Function Documentation	161
27.10.3.1	GetName	161
27.10.3.2	Print	161
27.10.3.3	Read	161
27.10.3.4	SetName	161
27.10.3.5	Size	161
27.10.3.6	Write	161
27.11	gdcm::ApplicationEntity Class Reference	161
27.11.1	Detailed Description	162
27.11.2	Member Function Documentation	162
27.11.2.1	IsValid	162
27.11.2.2	Print	162
27.11.2.3	SetBlob	162
27.11.2.4	Squeeze	162
27.11.3	Member Data Documentation	162
27.11.3.1	Internal	162
27.11.3.2	MaxLength	162
27.11.3.3	MaxNumberOfComponents	162
27.11.3.4	Padding	162
27.11.3.5	Separator	162
27.12	gdcm::network::AReleaseRPPDU Class Reference	163

27.12.1 Detailed Description	164
27.12.2 Constructor & Destructor Documentation	164
27.12.2.1 AReleaseRPPDU	164
27.12.3 Member Function Documentation	164
27.12.3.1 IsLastFragment	164
27.12.3.2 Print	164
27.12.3.3 Read	164
27.12.3.4 Size	164
27.12.3.5 Write	164
27.13gdcmm::network::AReleaseRQPDU Class Reference	164
27.13.1 Detailed Description	165
27.13.2 Constructor & Destructor Documentation	165
27.13.2.1 AReleaseRQPDU	165
27.13.3 Member Function Documentation	166
27.13.3.1 IsLastFragment	166
27.13.3.2 Print	166
27.13.3.3 Read	166
27.13.3.4 Size	166
27.13.3.5 Write	166
27.14gdcmm::network::ARTIMTimer Class Reference	166
27.14.1 Detailed Description	166
27.14.2 Constructor & Destructor Documentation	167
27.14.2.1 ARTIMTimer	167
27.14.3 Member Function Documentation	167
27.14.3.1 GetElapsedTime	167
27.14.3.2 GetHasExpired	167
27.14.3.3 GetTimeout	167
27.14.3.4 SetTimeout	167
27.14.3.5 Start	167
27.14.3.6 Stop	167
27.15gdcmm::ASN1 Class Reference	167
27.15.1 Detailed Description	168
27.15.2 Constructor & Destructor Documentation	168
27.15.2.1 ASN1	168
27.15.2.2 ~ASN1	168
27.15.3 Member Function Documentation	168
27.15.3.1 ParseDump	168

27.15.3.2 ParseDumpFile . . . . .	168
27.15.3.3 TestPBKDF2 . . . . .	168
27.16gdcmm::network::AsynchronousOperationsWindowSub Class Reference . . . . .	168
27.16.1 Detailed Description . . . . .	168
27.16.2 Constructor & Destructor Documentation . . . . .	168
27.16.2.1 AsynchronousOperationsWindowSub . . . . .	168
27.16.3 Member Function Documentation . . . . .	168
27.16.3.1 Print . . . . .	169
27.16.3.2 Read . . . . .	169
27.16.3.3 Size . . . . .	169
27.16.3.4 Write . . . . .	169
27.17gdcmm::Attribute< Group, Element, TVR, TVM > Class Template Reference . . . . .	169
27.17.1 Detailed Description . . . . .	170
27.17.2 Member Typedef Documentation . . . . .	171
27.17.2.1 ArrayType . . . . .	171
27.17.3 Member Enumeration Documentation . . . . .	171
27.17.3.1 anonymous enum . . . . .	171
27.17.4 Member Function Documentation . . . . .	171
27.17.4.1 GDCM_STATIC_ASSERT . . . . .	171
27.17.4.2 GDCM_STATIC_ASSERT . . . . .	171
27.17.4.3 GDCM_STATIC_ASSERT . . . . .	171
27.17.4.4 GetAsDataElement . . . . .	171
27.17.4.5 GetDictVM . . . . .	171
27.17.4.6 GetDictVR . . . . .	172
27.17.4.7 GetNumberOfValues . . . . .	172
27.17.4.8 GetTag . . . . .	172
27.17.4.9 GetValue . . . . .	172
27.17.4.10GetValue . . . . .	172
27.17.4.11GetValues . . . . .	173
27.17.4.12GetVM . . . . .	173
27.17.4.13GetVR . . . . .	173
27.17.4.14operator"!=" . . . . .	173
27.17.4.15operator< . . . . .	173
27.17.4.16operator== . . . . .	173
27.17.4.17operator[] . . . . .	174
27.17.4.18operator[] . . . . .	174
27.17.4.19Print . . . . .	174



27.17.4.20Set	174
27.17.4.21SetByteValue	174
27.17.4.22SetByteValueNoSwap	174
27.17.4.23SetFromDataElement	175
27.17.4.24SetFromDataSet	175
27.17.4.25SetValue	175
27.17.4.26SetValues	175
27.17.5 Member Data Documentation	175
27.17.5.1 Internal	175
27.18gdcmm::Attribute< Group, Element, TVR, VM::VM1 > Class Template Reference	176
27.18.1 Member Typedef Documentation	177
27.18.1.1 ArrayType	177
27.18.2 Member Enumeration Documentation	177
27.18.2.1 anonymous enum	177
27.18.3 Member Function Documentation	177
27.18.3.1 GDCM_STATIC_ASSERT	177
27.18.3.2 GDCM_STATIC_ASSERT	177
27.18.3.3 GDCM_STATIC_ASSERT	178
27.18.3.4 GDCM_STATIC_ASSERT	178
27.18.3.5 GetAsDataElement	178
27.18.3.6 GetDictVM	178
27.18.3.7 GetDictVR	178
27.18.3.8 GetNumberOfValues	178
27.18.3.9 GetTag	178
27.18.3.10GetValue	178
27.18.3.11GetValue	178
27.18.3.12GetValues	178
27.18.3.13GetVM	178
27.18.3.14GetVR	179
27.18.3.15operator!=	179
27.18.3.16operator<	179
27.18.3.17operator==	179
27.18.3.18Print	179
27.18.3.19Set	179
27.18.3.20SetByteValue	179
27.18.3.21SetByteValueNoSwap	179
27.18.3.22SetFromDataElement	180

27.18.3.23SetFromDataSet . . . . .	180
27.18.3.24SetValue . . . . .	180
27.18.4 Member Data Documentation . . . . .	180
27.18.4.1 Internal . . . . .	180
27.19gdcmm::Attribute< Group, Element, TVR, VM::VM1_3 > Class Template Reference . . . . .	180
27.19.1 Member Function Documentation . . . . .	181
27.19.1.1 GetVM . . . . .	181
27.20gdcmm::Attribute< Group, Element, TVR, VM::VM1_8 > Class Template Reference . . . . .	181
27.20.1 Member Function Documentation . . . . .	182
27.20.1.1 GetVM . . . . .	182
27.21gdcmm::Attribute< Group, Element, TVR, VM::VM1_n > Class Template Reference . . . . .	183
27.21.1 Member Typedef Documentation . . . . .	184
27.21.1.1 ArrayType . . . . .	184
27.21.2 Constructor & Destructor Documentation . . . . .	184
27.21.2.1 Attribute . . . . .	184
27.21.2.2 ~Attribute . . . . .	184
27.21.3 Member Function Documentation . . . . .	184
27.21.3.1 GDCM_STATIC_ASSERT . . . . .	184
27.21.3.2 GDCM_STATIC_ASSERT . . . . .	184
27.21.3.3 GDCM_STATIC_ASSERT . . . . .	184
27.21.3.4 GetAsDataElement . . . . .	184
27.21.3.5 GetDictVM . . . . .	185
27.21.3.6 GetDictVR . . . . .	185
27.21.3.7 GetNumberOfValues . . . . .	185
27.21.3.8 GetTag . . . . .	185
27.21.3.9 GetValue . . . . .	185
27.21.3.10GetValue . . . . .	185
27.21.3.11GetValues . . . . .	185
27.21.3.12GetVM . . . . .	185
27.21.3.13GetVR . . . . .	185
27.21.3.14operator[] . . . . .	185
27.21.3.15operator[] . . . . .	185
27.21.3.16Print . . . . .	186
27.21.3.17Set . . . . .	186
27.21.3.18SetByteValue . . . . .	186
27.21.3.19SetFromDataElement . . . . .	186
27.21.3.20SetFromDataSet . . . . .	186

27.21.3.21SetNumberOfValues . . . . .	186
27.21.3.22SetValue . . . . .	186
27.21.3.23SetValue . . . . .	186
27.21.3.24SetValues . . . . .	187
27.22gdcmm::Attribute< Group, Element, TVR, VM::VM2_n > Class Template Reference . . . . .	187
27.22.1 Member Function Documentation . . . . .	188
27.22.1.1 GetVM . . . . .	188
27.23gdcmm::Attribute< Group, Element, TVR, VM::VM2_n > Class Template Reference . . . . .	188
27.23.1 Member Function Documentation . . . . .	189
27.23.1.1 GetVM . . . . .	190
27.24gdcmm::Attribute< Group, Element, TVR, VM::VM3_3n > Class Template Reference . . . . .	190
27.24.1 Member Function Documentation . . . . .	191
27.24.1.1 GetVM . . . . .	191
27.25gdcmm::Attribute< Group, Element, TVR, VM::VM3_n > Class Template Reference . . . . .	191
27.25.1 Member Function Documentation . . . . .	192
27.25.1.1 GetVM . . . . .	193
27.26gdcmm::AudioCodec Class Reference . . . . .	193
27.26.1 Detailed Description . . . . .	194
27.26.2 Constructor & Destructor Documentation . . . . .	194
27.26.2.1 AudioCodec . . . . .	194
27.26.2.2 ~AudioCodec . . . . .	194
27.26.3 Member Function Documentation . . . . .	194
27.26.3.1 CanCode . . . . .	194
27.26.3.2 CanDecode . . . . .	195
27.26.3.3 Decode . . . . .	195
27.27gdcmm::Base64 Class Reference . . . . .	195
27.27.1 Detailed Description . . . . .	195
27.27.2 Member Function Documentation . . . . .	195
27.27.2.1 Decode . . . . .	195
27.27.2.2 Encode . . . . .	196
27.27.2.3 GetDecodeLength . . . . .	196
27.27.2.4 GetEncodeLength . . . . .	196
27.28gdcmm::network::BaseCompositeMessage Class Reference . . . . .	196
27.28.1 Detailed Description . . . . .	197
27.28.2 Constructor & Destructor Documentation . . . . .	198
27.28.2.1 ~BaseCompositeMessage . . . . .	198
27.28.3 Member Function Documentation . . . . .	198

27.28.3.1 ConstructPDV . . . . .	198
27.29gdcmm::network::BasePDU Class Reference . . . . .	198
27.29.1 Detailed Description . . . . .	199
27.29.2 Constructor & Destructor Documentation . . . . .	200
27.29.2.1 ~BasePDU . . . . .	200
27.29.3 Member Function Documentation . . . . .	200
27.29.3.1 IsLastFragment . . . . .	200
27.29.3.2 Print . . . . .	200
27.29.3.3 Read . . . . .	200
27.29.3.4 Size . . . . .	200
27.29.3.5 Write . . . . .	200
27.30gdcmm::BaseRootQuery Class Reference . . . . .	200
27.30.1 Detailed Description . . . . .	202
27.30.2 Constructor & Destructor Documentation . . . . .	202
27.30.2.1 BaseRootQuery . . . . .	202
27.30.2.2 ~BaseRootQuery . . . . .	202
27.30.3 Member Function Documentation . . . . .	202
27.30.3.1 AddQueryDataSet . . . . .	203
27.30.3.2 Construct . . . . .	203
27.30.3.3 GetAbstractSyntaxUID . . . . .	203
27.30.3.4 GetQueryDataSet . . . . .	203
27.30.3.5 GetQueryDataSet . . . . .	203
27.30.3.6 GetQueryLevelFromQueryRoot . . . . .	203
27.30.3.7 GetQueryLevelFromString . . . . .	203
27.30.3.8 GetQueryLevelString . . . . .	203
27.30.3.9 GetTagListByLevel . . . . .	203
27.30.3.10InitializeDataSet . . . . .	203
27.30.3.11Print . . . . .	203
27.30.3.12SetSearchParameter . . . . .	203
27.30.3.13SetSearchParameter . . . . .	203
27.30.3.14SetSearchParameter . . . . .	203
27.30.3.15ValidateQuery . . . . .	204
27.30.3.16WriteHelpFile . . . . .	204
27.30.3.17WriteQuery . . . . .	204
27.30.4 Friends And Related Function Documentation . . . . .	204
27.30.4.1 QueryFactory . . . . .	204
27.30.5 Member Data Documentation . . . . .	204

27.30.5.1 mDataSet . . . . .	204
27.30.5.2 mHelpDescription . . . . .	204
27.30.5.3 mImage . . . . .	204
27.30.5.4 mPatient . . . . .	204
27.30.5.5 mRootType . . . . .	204
27.30.5.6 mSeries . . . . .	204
27.30.5.7 mStudy . . . . .	204
27.31gdcm::SegmentHelper::BasicCodedEntry Struct Reference . . . . .	204
27.31.1 Detailed Description . . . . .	206
27.31.2 Constructor & Destructor Documentation . . . . .	206
27.31.2.1 BasicCodedEntry . . . . .	206
27.31.2.2 BasicCodedEntry . . . . .	206
27.31.2.3 BasicCodedEntry . . . . .	206
27.31.3 Member Function Documentation . . . . .	206
27.31.3.1 IsEmpty . . . . .	206
27.31.4 Member Data Documentation . . . . .	206
27.31.4.1 CM . . . . .	206
27.31.4.2 CSD . . . . .	206
27.31.4.3 CSV . . . . .	206
27.31.4.4 CV . . . . .	207
27.32gdcm::BasicOffsetTable Class Reference . . . . .	207
27.32.1 Detailed Description . . . . .	208
27.32.2 Constructor & Destructor Documentation . . . . .	208
27.32.2.1 BasicOffsetTable . . . . .	208
27.32.3 Member Function Documentation . . . . .	208
27.32.3.1 Read . . . . .	209
27.32.4 Friends And Related Function Documentation . . . . .	209
27.32.4.1 operator<< . . . . .	209
27.33gdcm::Bitmap Class Reference . . . . .	209
27.33.1 Detailed Description . . . . .	212
27.33.2 Member Typedef Documentation . . . . .	212
27.33.2.1 LUTPtr . . . . .	212
27.33.3 Constructor & Destructor Documentation . . . . .	212
27.33.3.1 Bitmap . . . . .	212
27.33.3.2 ~Bitmap . . . . .	212
27.33.4 Member Function Documentation . . . . .	212
27.33.4.1 AreOverlaysInPixelData . . . . .	212

27.33.4.2 Clear	212
27.33.4.3 ComputeLossyFlag	212
27.33.4.4 GetBuffer	212
27.33.4.5 GetBuffer2	212
27.33.4.6 GetBufferLength	212
27.33.4.7 GetColumns	213
27.33.4.8 GetDataElement	213
27.33.4.9 GetDataElement	213
27.33.4.10 GetDimension	213
27.33.4.11 GetDimensions	213
27.33.4.12 GetLUT	213
27.33.4.13 GetLUT	213
27.33.4.14 GetNeedByteSwap	213
27.33.4.15 GetNumberOfDimensions	213
27.33.4.16 GetPhotometricInterpretation	213
27.33.4.17 GetPixelFormat	214
27.33.4.18 GetPixelFormat	214
27.33.4.19 GetPlanarConfiguration	214
27.33.4.20 GetRows	214
27.33.4.21 GetTransferSyntax	214
27.33.4.22 IsEmpty	214
27.33.4.23 IsLossy	214
27.33.4.24 IsTransferSyntaxCompatible	214
27.33.4.25 Print	214
27.33.4.26 SetColumns	214
27.33.4.27 SetDataElement	214
27.33.4.28 SetDimension	215
27.33.4.29 SetDimensions	215
27.33.4.30 SetLossyFlag	215
27.33.4.31 SetLUT	215
27.33.4.32 SetNeedByteSwap	215
27.33.4.33 SetNumberOfDimensions	215
27.33.4.34 SetPhotometricInterpretation	215
27.33.4.35 SetPixelFormat	215
27.33.4.36 SetPlanarConfiguration	215
27.33.4.37 SetRows	216
27.33.4.38 SetTransferSyntax	216

27.33.4.39TryJPEG2000Codec . . . . .	216
27.33.4.40TryJPEG2000Codec2 . . . . .	216
27.33.4.41TryJPEGCodec . . . . .	216
27.33.4.42TryJPEGCodec2 . . . . .	216
27.33.4.43TryJPEGLSCodec . . . . .	216
27.33.4.44TryKAKADUCodec . . . . .	216
27.33.4.45TryPVRGCodec . . . . .	216
27.33.4.46TryRAWCodec . . . . .	216
27.33.4.47TryRLECodec . . . . .	216
27.33.5 Friends And Related Function Documentation . . . . .	216
27.33.5.1 ImageChangeTransferSyntax . . . . .	216
27.33.5.2 PixmapReader . . . . .	216
27.33.6 Member Data Documentation . . . . .	216
27.33.6.1 Dimensions . . . . .	216
27.33.6.2 LossyFlag . . . . .	216
27.33.6.3 LUT . . . . .	216
27.33.6.4 NeedByteSwap . . . . .	216
27.33.6.5 NumberOfDimensions . . . . .	216
27.33.6.6 PF . . . . .	216
27.33.6.7 PI . . . . .	216
27.33.6.8 PixelData . . . . .	216
27.33.6.9 PlanarConfiguration . . . . .	217
27.33.6.10TS . . . . .	217
27.34gdcm::BitmapToBitmapFilter Class Reference . . . . .	217
27.34.1 Detailed Description . . . . .	218
27.34.2 Constructor & Destructor Documentation . . . . .	218
27.34.2.1 BitmapToBitmapFilter . . . . .	218
27.34.2.2 ~BitmapToBitmapFilter . . . . .	218
27.34.3 Member Function Documentation . . . . .	218
27.34.3.1 GetOutput . . . . .	218
27.34.3.2 GetOutputAsBitmap . . . . .	218
27.34.3.3 SetInput . . . . .	218
27.34.4 Member Data Documentation . . . . .	218
27.34.4.1 Input . . . . .	218
27.34.4.2 Output . . . . .	218
27.35gdcm::BoxRegion Class Reference . . . . .	219
27.35.1 Detailed Description . . . . .	220

27.35.2 Constructor & Destructor Documentation . . . . .	220
27.35.2.1 BoxRegion . . . . .	220
27.35.2.2 ~BoxRegion . . . . .	220
27.35.2.3 BoxRegion . . . . .	220
27.35.3 Member Function Documentation . . . . .	220
27.35.3.1 Area . . . . .	220
27.35.3.2 BoundingBox . . . . .	221
27.35.3.3 Clone . . . . .	221
27.35.3.4 ComputeBoundingBox . . . . .	221
27.35.3.5 Empty . . . . .	221
27.35.3.6 GetXMax . . . . .	221
27.35.3.7 GetXMin . . . . .	221
27.35.3.8 GetYMax . . . . .	221
27.35.3.9 GetYMin . . . . .	221
27.35.3.10GetZMax . . . . .	221
27.35.3.11GetZMin . . . . .	221
27.35.3.12IsValid . . . . .	221
27.35.3.13operator= . . . . .	221
27.35.3.14Print . . . . .	221
27.35.3.15SetDomain . . . . .	222
27.36gdcmm::ByteBuffer Class Reference . . . . .	222
27.36.1 Detailed Description . . . . .	222
27.36.2 Constructor & Destructor Documentation . . . . .	222
27.36.2.1 ByteBuffer . . . . .	222
27.36.3 Member Function Documentation . . . . .	222
27.36.3.1 Get . . . . .	222
27.36.3.2 GetStart . . . . .	222
27.36.3.3 ShiftEnd . . . . .	222
27.36.3.4 UpdatePosition . . . . .	222
27.37gdcmm::ByteSwap< T > Class Template Reference . . . . .	223
27.37.1 Detailed Description . . . . .	223
27.37.2 Member Function Documentation . . . . .	223
27.37.2.1 Swap . . . . .	223
27.37.2.2 SwapFromSwapCodeIntoSystem . . . . .	223
27.37.2.3 SwapRange . . . . .	223
27.37.2.4 SwapRangeFromSwapCodeIntoSystem . . . . .	223
27.37.2.5 SystemIsBigEndian . . . . .	223



27.37.2.6 SystemIsLittleEndian . . . . .	224
27.38gdcmm::ByteSwapFilter Class Reference . . . . .	224
27.38.1 Detailed Description . . . . .	224
27.38.2 Constructor & Destructor Documentation . . . . .	224
27.38.2.1 ByteSwapFilter . . . . .	224
27.38.2.2 ~ByteSwapFilter . . . . .	224
27.38.3 Member Function Documentation . . . . .	224
27.38.3.1 ByteSwap . . . . .	224
27.38.3.2 SetByteSwapTag . . . . .	224
27.39gdcmm::ByteValue Class Reference . . . . .	224
27.39.1 Detailed Description . . . . .	226
27.39.2 Constructor & Destructor Documentation . . . . .	226
27.39.2.1 ByteValue . . . . .	226
27.39.2.2 ByteValue . . . . .	227
27.39.2.3 ~ByteValue . . . . .	227
27.39.3 Member Function Documentation . . . . .	227
27.39.3.1 Append . . . . .	227
27.39.3.2 Clear . . . . .	227
27.39.3.3 ComputeLength . . . . .	227
27.39.3.4 Fill . . . . .	227
27.39.3.5 GetBuffer . . . . .	227
27.39.3.6 GetLength . . . . .	227
27.39.3.7 GetPointer . . . . .	228
27.39.3.8 IsEmpty . . . . .	228
27.39.3.9 IsPrintable . . . . .	228
27.39.3.10operator const std::vector< char > & . . . . .	228
27.39.3.11operator= . . . . .	228
27.39.3.12operator== . . . . .	228
27.39.3.13operator== . . . . .	228
27.39.3.14Print . . . . .	228
27.39.3.15PrintASCII . . . . .	228
27.39.3.16PrintASCIIXML . . . . .	228
27.39.3.17PrintGroupLength . . . . .	228
27.39.3.18PrintHex . . . . .	228
27.39.3.19PrintHexXML . . . . .	228
27.39.3.20PrintPNXML . . . . .	228
27.39.3.21Read . . . . .	229

27.39.3.22Read . . . . .	229
27.39.3.23SetLength . . . . .	229
27.39.3.24SetLengthOnly . . . . .	229
27.39.3.25Write . . . . .	229
27.39.3.26Write . . . . .	229
27.39.3.27WriteBuffer . . . . .	229
27.40gdcm::CAPICryptoFactory Class Reference . . . . .	229
27.40.1 Constructor & Destructor Documentation . . . . .	230
27.40.1.1 CAPICryptoFactory . . . . .	230
27.40.2 Member Function Documentation . . . . .	230
27.40.2.1 CreateCMSProvider . . . . .	230
27.41gdcm::CAPICryptographicMessageSyntax Class Reference . . . . .	230
27.41.1 Constructor & Destructor Documentation . . . . .	232
27.41.1.1 CAPICryptographicMessageSyntax . . . . .	232
27.41.1.2 ~CAPICryptographicMessageSyntax . . . . .	232
27.41.2 Member Function Documentation . . . . .	232
27.41.2.1 Decrypt . . . . .	232
27.41.2.2 Encrypt . . . . .	232
27.41.2.3 GetCipherType . . . . .	232
27.41.2.4 GetInitialized . . . . .	232
27.41.2.5 ParseCertificateFile . . . . .	232
27.41.2.6 ParseKeyFile . . . . .	232
27.41.2.7 SetCipherType . . . . .	232
27.41.2.8 SetPassword . . . . .	232
27.42gdcm::network::CEchoRQ Class Reference . . . . .	233
27.42.1 Detailed Description . . . . .	234
27.42.2 Member Function Documentation . . . . .	234
27.42.2.1 ConstructPDV . . . . .	234
27.42.3 Member Data Documentation . . . . .	234
27.42.3.1 AffectedSOPClassUID . . . . .	234
27.42.3.2 MessageID . . . . .	234
27.43gdcm::network::CEchoRSP Class Reference . . . . .	234
27.43.1 Detailed Description . . . . .	235
27.43.2 Member Function Documentation . . . . .	235
27.43.2.1 ConstructPDVByDataSet . . . . .	235
27.44gdcm::network::CFind Class Reference . . . . .	235
27.44.1 Detailed Description . . . . .	235

27.45gdcmm::network::CFindCancelRQ Class Reference . . . . .	236
27.45.1 Detailed Description . . . . .	236
27.45.2 Member Function Documentation . . . . .	236
27.45.2.1 ConstructPDVByDataSet . . . . .	237
27.46gdcmm::network::CFindRQ Class Reference . . . . .	237
27.46.1 Detailed Description . . . . .	238
27.46.2 Member Function Documentation . . . . .	238
27.46.2.1 ConstructPDV . . . . .	238
27.47gdcmm::network::CFindRSP Class Reference . . . . .	238
27.47.1 Detailed Description . . . . .	239
27.47.2 Member Function Documentation . . . . .	239
27.47.2.1 ConstructPDVByDataSet . . . . .	239
27.48gdcmm::network::CMoveCancelRq Class Reference . . . . .	239
27.48.1 Member Function Documentation . . . . .	240
27.48.1.1 ConstructPDVByDataSet . . . . .	240
27.49gdcmm::network::CMoveRQ Class Reference . . . . .	241
27.49.1 Detailed Description . . . . .	241
27.49.2 Member Function Documentation . . . . .	241
27.49.2.1 ConstructPDV . . . . .	242
27.50gdcmm::network::CMoveRSP Class Reference . . . . .	242
27.50.1 Detailed Description . . . . .	243
27.50.2 Member Function Documentation . . . . .	243
27.50.2.1 ConstructPDVByDataSet . . . . .	243
27.51gdcmm::Codec Class Reference . . . . .	243
27.51.1 Detailed Description . . . . .	244
27.52gdcmm::Coder Class Reference . . . . .	244
27.52.1 Detailed Description . . . . .	245
27.52.2 Constructor & Destructor Documentation . . . . .	245
27.52.2.1 ~Coder . . . . .	245
27.52.3 Member Function Documentation . . . . .	245
27.52.3.1 CanCode . . . . .	245
27.52.3.2 Code . . . . .	245
27.52.3.3 InternalCode . . . . .	245
27.53gdcmm::CodeString Class Reference . . . . .	245
27.53.1 Detailed Description . . . . .	246
27.53.2 Member Typedef Documentation . . . . .	247
27.53.2.1 const_iterator . . . . .	247

27.53.2.2 const_reference . . . . .	247
27.53.2.3 const_reverse_iterator . . . . .	247
27.53.2.4 difference_type . . . . .	247
27.53.2.5 iterator . . . . .	247
27.53.2.6 pointer . . . . .	247
27.53.2.7 reference . . . . .	247
27.53.2.8 reverse_iterator . . . . .	247
27.53.2.9 size_type . . . . .	247
27.53.2.10value_type . . . . .	247
27.53.3 Constructor & Destructor Documentation . . . . .	247
27.53.3.1 CodeString . . . . .	247
27.53.3.2 CodeString . . . . .	247
27.53.3.3 CodeString . . . . .	247
27.53.3.4 CodeString . . . . .	247
27.53.4 Member Function Documentation . . . . .	247
27.53.4.1 GetAsString . . . . .	248
27.53.4.2 IsValid . . . . .	248
27.53.4.3 Size . . . . .	248
27.53.4.4 TrimInternal . . . . .	248
27.53.5 Friends And Related Function Documentation . . . . .	248
27.53.5.1 operator"!=" . . . . .	248
27.53.5.2 operator<< . . . . .	248
27.53.5.3 operator== . . . . .	248
27.54gdcmm::Command Class Reference . . . . .	248
27.54.1 Detailed Description . . . . .	250
27.54.2 Constructor & Destructor Documentation . . . . .	250
27.54.2.1 Command . . . . .	250
27.54.2.2 ~Command . . . . .	250
27.54.3 Member Function Documentation . . . . .	250
27.54.3.1 Execute . . . . .	250
27.54.3.2 Execute . . . . .	250
27.55gdcmm::CommandDataSet Class Reference . . . . .	250
27.55.1 Detailed Description . . . . .	252
27.55.2 Constructor & Destructor Documentation . . . . .	252
27.55.2.1 CommandDataSet . . . . .	252
27.55.2.2 ~CommandDataSet . . . . .	252
27.55.3 Member Function Documentation . . . . .	252

27.55.3.1 Insert . . . . .	252
27.55.3.2 Read . . . . .	252
27.55.3.3 Replace . . . . .	252
27.55.3.4 Write . . . . .	252
27.55.4 Friends And Related Function Documentation . . . . .	252
27.55.4.1 operator<< . . . . .	252
27.56gdcmm::network::CompositeMessageFactory Class Reference . . . . .	252
27.56.1 Detailed Description . . . . .	253
27.56.2 Member Function Documentation . . . . .	253
27.56.2.1 ConstructCEchoRQ . . . . .	253
27.56.2.2 ConstructCFindRQ . . . . .	253
27.56.2.3 ConstructCMoveRQ . . . . .	253
27.56.2.4 ConstructCStoreRQ . . . . .	253
27.56.2.5 ConstructCStoreRSP . . . . .	253
27.57gdcmm::CompositeNetworkFunctions Class Reference . . . . .	253
27.57.1 Detailed Description . . . . .	254
27.57.2 Member Typedef Documentation . . . . .	254
27.57.2.1 KeyValuePairArrayType . . . . .	254
27.57.2.2 KeyValuePairType . . . . .	254
27.57.3 Member Function Documentation . . . . .	254
27.57.3.1 CEcho . . . . .	255
27.57.3.2 CFind . . . . .	256
27.57.3.3 CMove . . . . .	256
27.57.3.4 ConstructQuery . . . . .	257
27.57.3.5 ConstructQuery . . . . .	257
27.57.3.6 CStore . . . . .	257
27.58gdcmm::ConstCharWrapper Class Reference . . . . .	258
27.58.1 Detailed Description . . . . .	258
27.58.2 Constructor & Destructor Documentation . . . . .	258
27.58.2.1 ConstCharWrapper . . . . .	258
27.58.3 Member Function Documentation . . . . .	258
27.58.3.1 operator const char * . . . . .	258
27.59gdcmm::CP246ExplicitDataElement Class Reference . . . . .	258
27.59.1 Detailed Description . . . . .	259
27.59.2 Member Function Documentation . . . . .	259
27.59.2.1 GetLength . . . . .	259
27.59.2.2 Read . . . . .	259

27.59.2.3 ReadPreValue . . . . .	260
27.59.2.4 ReadValue . . . . .	260
27.59.2.5 ReadWithLength . . . . .	260
27.60gdcmm::CryptoFactory Class Reference . . . . .	260
27.60.1 Detailed Description . . . . .	261
27.60.2 Member Enumeration Documentation . . . . .	261
27.60.2.1 CryptoLib . . . . .	261
27.60.3 Constructor & Destructor Documentation . . . . .	261
27.60.3.1 CryptoFactory . . . . .	261
27.60.3.2 CryptoFactory . . . . .	261
27.60.3.3 ~CryptoFactory . . . . .	261
27.60.4 Member Function Documentation . . . . .	261
27.60.4.1 CreateCMSProvider . . . . .	261
27.60.4.2 GetFactoryInstance . . . . .	261
27.61gdcmm::CryptographicMessageSyntax Class Reference . . . . .	262
27.61.1 Member Enumeration Documentation . . . . .	262
27.61.1.1 CipherTypes . . . . .	262
27.61.2 Constructor & Destructor Documentation . . . . .	263
27.61.2.1 CryptographicMessageSyntax . . . . .	263
27.61.2.2 ~CryptographicMessageSyntax . . . . .	263
27.61.3 Member Function Documentation . . . . .	263
27.61.3.1 Decrypt . . . . .	263
27.61.3.2 Encrypt . . . . .	263
27.61.3.3 GetCipherType . . . . .	263
27.61.3.4 ParseCertificateFile . . . . .	263
27.61.3.5 ParseKeyFile . . . . .	263
27.61.3.6 SetCipherType . . . . .	263
27.61.3.7 SetPassword . . . . .	264
27.62gdcmm::CSAElement Class Reference . . . . .	264
27.62.1 Detailed Description . . . . .	265
27.62.2 Member Typedef Documentation . . . . .	266
27.62.2.1 DataPtr . . . . .	266
27.62.3 Constructor & Destructor Documentation . . . . .	266
27.62.3.1 CSAElement . . . . .	266
27.62.3.2 CSAElement . . . . .	266
27.62.4 Member Function Documentation . . . . .	266
27.62.4.1 GetByteValue . . . . .	266

27.62.4.2 GetKey . . . . .	266
27.62.4.3 GetName . . . . .	266
27.62.4.4 GetNoOfItems . . . . .	266
27.62.4.5 GetSyngoDT . . . . .	266
27.62.4.6 GetValue . . . . .	266
27.62.4.7 GetValue . . . . .	267
27.62.4.8 GetVM . . . . .	267
27.62.4.9 GetVR . . . . .	267
27.62.4.10 IsEmpty . . . . .	267
27.62.4.11 operator< . . . . .	267
27.62.4.12 operator= . . . . .	267
27.62.4.13 operator== . . . . .	267
27.62.4.14 SetByteValue . . . . .	267
27.62.4.15 SetKey . . . . .	267
27.62.4.16 SetName . . . . .	267
27.62.4.17 SetNoOfItems . . . . .	267
27.62.4.18 SetSyngoDT . . . . .	267
27.62.4.19 SetValue . . . . .	267
27.62.4.20 SetVM . . . . .	267
27.62.4.21 SetVR . . . . .	267
27.62.5 Friends And Related Function Documentation . . . . .	268
27.62.5.1 operator<< . . . . .	268
27.62.6 Member Data Documentation . . . . .	268
27.62.6.1 DataField . . . . .	268
27.62.6.2 KeyField . . . . .	268
27.62.6.3 NameField . . . . .	268
27.62.6.4 NoOfItemsField . . . . .	268
27.62.6.5 SyngoDTField . . . . .	268
27.62.6.6 ValueMultiplicityField . . . . .	268
27.62.6.7 VRField . . . . .	268
27.63gdcmm::CSAHeader Class Reference . . . . .	268
27.63.1 Detailed Description . . . . .	270
27.63.2 Member Enumeration Documentation . . . . .	270
27.63.2.1 CSAHeaderType . . . . .	270
27.63.3 Constructor & Destructor Documentation . . . . .	270
27.63.3.1 CSAHeader . . . . .	270
27.63.3.2 ~CSAHeader . . . . .	270

27.63.4 Member Function Documentation	270
27.63.4.1 FindCSAElementByName	271
27.63.4.2 GetCSADataInfo	271
27.63.4.3 GetCSAEnd	271
27.63.4.4 GetCSAElementByName	271
27.63.4.5 GetCSAImageHeaderInfoTag	271
27.63.4.6 GetCSASeriesHeaderInfoTag	271
27.63.4.7 GetDataSet	272
27.63.4.8 GetFormat	272
27.63.4.9 GetInterfile	272
27.63.4.10 LoadFromDataElement	272
27.63.4.11 Print	272
27.63.4.12 Read	272
27.63.4.13 Write	272
27.63.5 Friends And Related Function Documentation	272
27.63.5.1 operator<<	272
27.64 gdcmm::CSAHeaderDict Class Reference	272
27.64.1 Detailed Description	273
27.64.2 Member Typedef Documentation	273
27.64.2.1 ConstIterator	273
27.64.2.2 Iterator	273
27.64.2.3 MapCSAHeaderDictEntry	273
27.64.3 Constructor & Destructor Documentation	273
27.64.3.1 CSAHeaderDict	273
27.64.4 Member Function Documentation	273
27.64.4.1 AddCSAHeaderDictEntry	273
27.64.4.2 Begin	274
27.64.4.3 End	274
27.64.4.4 GetCSAHeaderDictEntry	274
27.64.4.5 IsEmpty	274
27.64.4.6 LoadDefault	274
27.64.5 Friends And Related Function Documentation	274
27.64.5.1 Dicts	274
27.64.5.2 operator<<	274
27.65 gdcmm::CSAHeaderDictEntry Class Reference	274
27.65.1 Detailed Description	275
27.65.2 Constructor & Destructor Documentation	275



27.65.2.1 CSAHeaderDictEntry . . . . .	275
27.65.3 Member Function Documentation . . . . .	275
27.65.3.1 GetDescription . . . . .	275
27.65.3.2 GetName . . . . .	275
27.65.3.3 GetVM . . . . .	275
27.65.3.4 GetVR . . . . .	275
27.65.3.5 operator< . . . . .	275
27.65.3.6 SetDescription . . . . .	276
27.65.3.7 SetName . . . . .	276
27.65.3.8 SetVM . . . . .	276
27.65.3.9 SetVR . . . . .	276
27.65.4 Friends And Related Function Documentation . . . . .	276
27.65.4.1 operator<< . . . . .	276
27.66gdcm::CSAHeaderDictException Class Reference . . . . .	276
27.67gdcm::network::CStoreRQ Class Reference . . . . .	277
27.67.1 Detailed Description . . . . .	278
27.67.2 Member Function Documentation . . . . .	278
27.67.2.1 ConstructPDV . . . . .	278
27.68gdcm::network::CStoreRSP Class Reference . . . . .	278
27.68.1 Detailed Description . . . . .	279
27.68.2 Member Function Documentation . . . . .	279
27.68.2.1 ConstructPDV . . . . .	279
27.69gdcm::Curve Class Reference . . . . .	280
27.69.1 Detailed Description . . . . .	281
27.69.2 Constructor & Destructor Documentation . . . . .	281
27.69.2.1 Curve . . . . .	281
27.69.2.2 ~Curve . . . . .	281
27.69.2.3 Curve . . . . .	281
27.69.3 Member Function Documentation . . . . .	281
27.69.3.1 Decode . . . . .	281
27.69.3.2 GetAsPoints . . . . .	281
27.69.3.3 GetCurveDataDescriptor . . . . .	282
27.69.3.4 GetDataValueRepresentation . . . . .	282
27.69.3.5 GetDimensions . . . . .	282
27.69.3.6 GetGroup . . . . .	282
27.69.3.7 GetNumberOfCurves . . . . .	282
27.69.3.8 GetNumberOfPoints . . . . .	282

27.69.3.9 GetTypeOfData . . . . .	282
27.69.3.10GetTypeOfDataDescription . . . . .	282
27.69.3.11IsEmpty . . . . .	282
27.69.3.12Print . . . . .	282
27.69.3.13SetCoordinateStartValue . . . . .	282
27.69.3.14SetCoordinateStepValue . . . . .	282
27.69.3.15SetCurve . . . . .	282
27.69.3.16SetCurveDataDescriptor . . . . .	282
27.69.3.17SetCurveDescription . . . . .	282
27.69.3.18SetDataValueRepresentation . . . . .	282
27.69.3.19SetDimensions . . . . .	282
27.69.3.20SetGroup . . . . .	282
27.69.3.21SetNumberOfPoints . . . . .	282
27.69.3.22SetTypeOfData . . . . .	282
27.69.3.23Update . . . . .	282
27.70gdcmm::DataElement Class Reference . . . . .	282
27.70.1 Detailed Description . . . . .	285
27.70.2 Member Typedef Documentation . . . . .	285
27.70.2.1 ValuePtr . . . . .	285
27.70.3 Constructor & Destructor Documentation . . . . .	286
27.70.3.1 DataElement . . . . .	286
27.70.3.2 DataElement . . . . .	286
27.70.4 Member Function Documentation . . . . .	286
27.70.4.1 Clear . . . . .	286
27.70.4.2 Empty . . . . .	286
27.70.4.3 GetByteValue . . . . .	286
27.70.4.4 GetLength . . . . .	286
27.70.4.5 GetSequenceOfFragments . . . . .	286
27.70.4.6 GetSequenceOfFragments . . . . .	287
27.70.4.7 GetTag . . . . .	287
27.70.4.8 GetTag . . . . .	287
27.70.4.9 GetValue . . . . .	287
27.70.4.10GetValue . . . . .	287
27.70.4.11GetValueAsSQ . . . . .	287
27.70.4.12GetVL . . . . .	287
27.70.4.13GetVL . . . . .	288
27.70.4.14GetVR . . . . .	288

27.70.4.15	IsEmpty	288
27.70.4.16	IsUndefinedLength	288
27.70.4.17	operator<	288
27.70.4.18	operator=	288
27.70.4.19	operator==	288
27.70.4.20	Read	289
27.70.4.21	ReadOrSkip	289
27.70.4.22	ReadPreValue	289
27.70.4.23	ReadValue	289
27.70.4.24	ReadValueWithLength	289
27.70.4.25	ReadWithLength	289
27.70.4.26	SetByteValue	289
27.70.4.27	SetTag	289
27.70.4.28	SetValue	289
27.70.4.29	SetValueFieldLength	290
27.70.4.30	SetVL	290
27.70.4.31	SetVLToUndefined	290
27.70.4.32	SetVR	290
27.70.4.33	Write	291
27.70.5	Friends And Related Function Documentation	291
27.70.5.1	operator<<	291
27.70.6	Member Data Documentation	291
27.70.6.1	TagField	291
27.70.6.2	ValueField	291
27.70.6.3	ValueLengthField	291
27.70.6.4	VRField	291
27.71	gdcm::DataElementException Class Reference	291
27.72	gdcm::DataEvent Class Reference	292
27.72.1	Detailed Description	293
27.72.2	Member Typedef Documentation	293
27.72.2.1	Self	293
27.72.2.2	Superclass	293
27.72.3	Constructor & Destructor Documentation	293
27.72.3.1	DataEvent	294
27.72.3.2	~DataEvent	294
27.72.3.3	DataEvent	294
27.72.4	Member Function Documentation	294

27.72.4.1 CheckEvent . . . . .	294
27.72.4.2 GetData . . . . .	294
27.72.4.3 GetDataLength . . . . .	294
27.72.4.4 GetEventName . . . . .	294
27.72.4.5 MakeObject . . . . .	294
27.72.4.6 SetData . . . . .	294
27.73gdcm::DataSet Class Reference . . . . .	294
27.73.1 Detailed Description . . . . .	297
27.73.2 Member Typedef Documentation . . . . .	297
27.73.2.1 ConstIterator . . . . .	297
27.73.2.2 DataElementSet . . . . .	297
27.73.2.3 Iterator . . . . .	297
27.73.2.4 SizeType . . . . .	297
27.73.3 Member Function Documentation . . . . .	297
27.73.3.1 Begin . . . . .	297
27.73.3.2 Begin . . . . .	297
27.73.3.3 Clear . . . . .	297
27.73.3.4 ComputeDataElement . . . . .	298
27.73.3.5 ComputeGroupLength . . . . .	298
27.73.3.6 End . . . . .	298
27.73.3.7 End . . . . .	298
27.73.3.8 FindDataElement . . . . .	298
27.73.3.9 FindDataElement . . . . .	298
27.73.3.10FindNextDataElement . . . . .	298
27.73.3.11GetDataElement . . . . .	298
27.73.3.12GetDataElement . . . . .	299
27.73.3.13GetDEEnd . . . . .	299
27.73.3.14GetDES . . . . .	299
27.73.3.15GetDES . . . . .	299
27.73.3.16GetLength . . . . .	299
27.73.3.17GetMediaStorage . . . . .	299
27.73.3.18GetPrivateCreator . . . . .	299
27.73.3.19Insert . . . . .	299
27.73.3.20InsertDataElement . . . . .	299
27.73.3.21IsEmpty . . . . .	300
27.73.3.22operator() . . . . .	300
27.73.3.23operator= . . . . .	300

27.73.3.24operator[]	300
27.73.3.25Print	300
27.73.3.26Read	300
27.73.3.27ReadNested	300
27.73.3.28ReadSelectedPrivateTags	300
27.73.3.29ReadSelectedPrivateTagsWithLength	300
27.73.3.30ReadSelectedTags	300
27.73.3.31ReadSelectedTagsWithLength	300
27.73.3.32ReadUpToTag	300
27.73.3.33ReadUpToTagWithLength	300
27.73.3.34ReadWithLength	300
27.73.3.35Remove	300
27.73.3.36Replace	301
27.73.3.37ReplaceEmpty	301
27.73.3.38Size	301
27.73.3.39Write	301
27.73.4 Friends And Related Function Documentation	301
27.73.4.1 CSAHeader	301
27.73.4.2 operator<<	301
27.74gdcm::DataSetEvent Class Reference	301
27.74.1 Detailed Description	303
27.74.2 Member Typedef Documentation	303
27.74.2.1 Self	303
27.74.2.2 Superclass	303
27.74.3 Constructor & Destructor Documentation	303
27.74.3.1 DataSetEvent	303
27.74.3.2 ~DataSetEvent	303
27.74.3.3 DataSetEvent	303
27.74.4 Member Function Documentation	303
27.74.4.1 CheckEvent	303
27.74.4.2 GetDataSet	303
27.74.4.3 GetEventName	303
27.74.4.4 MakeObject	303
27.75gdcm::DataSetHelper Class Reference	304
27.75.1 Detailed Description	304
27.75.2 Member Function Documentation	304
27.75.2.1 ComputeVR	304

27.76gdcmm::Decoder Class Reference . . . . .	304
27.76.1 Detailed Description . . . . .	305
27.76.2 Constructor & Destructor Documentation . . . . .	305
27.76.2.1 ~Decoder . . . . .	305
27.76.3 Member Function Documentation . . . . .	305
27.76.3.1 CanDecode . . . . .	305
27.76.3.2 Decode . . . . .	305
27.76.3.3 DecodeByStreams . . . . .	305
27.77gdcmm::DefinedTerms Class Reference . . . . .	306
27.77.1 Detailed Description . . . . .	306
27.77.2 Constructor & Destructor Documentation . . . . .	306
27.77.2.1 DefinedTerms . . . . .	306
27.78gdcmm::Defs Class Reference . . . . .	306
27.78.1 Detailed Description . . . . .	307
27.78.2 Constructor & Destructor Documentation . . . . .	307
27.78.2.1 Defs . . . . .	307
27.78.2.2 ~Defs . . . . .	307
27.78.3 Member Function Documentation . . . . .	307
27.78.3.1 GetIODFromFile . . . . .	307
27.78.3.2 GetIODNameFromMediaStorage . . . . .	307
27.78.3.3 GetIODs . . . . .	308
27.78.3.4 GetIODs . . . . .	308
27.78.3.5 GetMacros . . . . .	308
27.78.3.6 GetMacros . . . . .	308
27.78.3.7 GetModules . . . . .	308
27.78.3.8 GetModules . . . . .	308
27.78.3.9 GetTypeFromTag . . . . .	308
27.78.3.10IsEmpty . . . . .	308
27.78.3.11LoadDefaults . . . . .	308
27.78.3.12LoadFromFile . . . . .	308
27.78.3.13Verify . . . . .	308
27.78.3.14Verify . . . . .	308
27.78.4 Friends And Related Function Documentation . . . . .	308
27.78.4.1 Global . . . . .	308
27.79gdcmm::DeltaEncodingCodec Class Reference . . . . .	309
27.79.1 Detailed Description . . . . .	310
27.79.2 Constructor & Destructor Documentation . . . . .	310

27.79.2.1 DeltaEncodingCodec	310
27.79.2.2 ~DeltaEncodingCodec	310
27.79.3 Member Function Documentation	310
27.79.3.1 CanDecode	310
27.79.3.2 Decode	310
27.79.3.3 Decode	310
27.80gdcm::DICOMDIR Class Reference	310
27.80.1 Detailed Description	310
27.80.2 Constructor & Destructor Documentation	311
27.80.2.1 DICOMDIR	311
27.80.2.2 DICOMDIR	311
27.81gdcm::DICOMDIRGenerator Class Reference	311
27.81.1 Detailed Description	312
27.81.2 Member Typedef Documentation	312
27.81.2.1 FilenamesType	312
27.81.2.2 FilenameType	312
27.81.3 Constructor & Destructor Documentation	312
27.81.3.1 DICOMDIRGenerator	312
27.81.3.2 ~DICOMDIRGenerator	312
27.81.4 Member Function Documentation	312
27.81.4.1 AddImageDirectoryRecord	312
27.81.4.2 AddPatientDirectoryRecord	312
27.81.4.3 AddSeriesDirectoryRecord	312
27.81.4.4 AddStudyDirectoryRecord	312
27.81.4.5 Generate	312
27.81.4.6 GetFile	313
27.81.4.7 GetScanner	313
27.81.4.8 SetDescriptor	313
27.81.4.9 SetFile	313
27.81.4.10SetFilenames	313
27.81.4.11SetRootDirectory	313
27.82gdcm::Dict Class Reference	313
27.82.1 Detailed Description	314
27.82.2 Member Typedef Documentation	314
27.82.2.1 ConstIterator	314
27.82.2.2 Iterator	314
27.82.2.3 MapDictEntry	314

27.82.3 Constructor & Destructor Documentation . . . . .	314
27.82.3.1 Dict . . . . .	314
27.82.4 Member Function Documentation . . . . .	314
27.82.4.1 AddDictEntry . . . . .	314
27.82.4.2 Begin . . . . .	314
27.82.4.3 End . . . . .	315
27.82.4.4 GetDictEntry . . . . .	315
27.82.4.5 GetDictEntryByKeyword . . . . .	315
27.82.4.6 GetDictEntryByName . . . . .	315
27.82.4.7 GetKeywordFromTag . . . . .	315
27.82.4.8 IsEmpty . . . . .	315
27.82.4.9 LoadDefault . . . . .	315
27.82.5 Friends And Related Function Documentation . . . . .	315
27.82.5.1 Dicts . . . . .	315
27.82.5.2 operator<< . . . . .	315
27.83gdcmm::DictConverter Class Reference . . . . .	316
27.83.1 Detailed Description . . . . .	316
27.83.2 Member Enumeration Documentation . . . . .	317
27.83.2.1 OutputTypes . . . . .	317
27.83.3 Constructor & Destructor Documentation . . . . .	317
27.83.3.1 DictConverter . . . . .	317
27.83.3.2 ~DictConverter . . . . .	317
27.83.4 Member Function Documentation . . . . .	317
27.83.4.1 AddGroupLength . . . . .	317
27.83.4.2 Convert . . . . .	317
27.83.4.3 ConvertToCXX . . . . .	317
27.83.4.4 ConvertToXML . . . . .	317
27.83.4.5 GetDictName . . . . .	317
27.83.4.6 GetInputFilename . . . . .	317
27.83.4.7 GetOutputFilename . . . . .	317
27.83.4.8 GetOutputType . . . . .	317
27.83.4.9 Readuint16 . . . . .	317
27.83.4.10ReadVM . . . . .	317
27.83.4.11ReadVR . . . . .	317
27.83.4.12SetDictName . . . . .	317
27.83.4.13SetInputFileName . . . . .	317
27.83.4.14SetOutputFileName . . . . .	317



27.83.4.15SetOutputType . . . . .	317
27.83.4.16WriteFooter . . . . .	318
27.83.4.17WriteHeader . . . . .	318
27.84gdcmm::DictEntry Class Reference . . . . .	318
27.84.1 Detailed Description . . . . .	319
27.84.2 Constructor & Destructor Documentation . . . . .	319
27.84.2.1 DictEntry . . . . .	319
27.84.3 Member Function Documentation . . . . .	319
27.84.3.1 GetKeyword . . . . .	319
27.84.3.2 GetName . . . . .	319
27.84.3.3 GetRetired . . . . .	319
27.84.3.4 GetVM . . . . .	319
27.84.3.5 GetVR . . . . .	320
27.84.3.6 IsUnique . . . . .	320
27.84.3.7 SetElementXX . . . . .	320
27.84.3.8 SetGroupXX . . . . .	320
27.84.3.9 SetKeyword . . . . .	320
27.84.3.10SetName . . . . .	320
27.84.3.11SetRetired . . . . .	320
27.84.3.12SetVM . . . . .	320
27.84.3.13SetVR . . . . .	320
27.84.4 Friends And Related Function Documentation . . . . .	320
27.84.4.1 Dict . . . . .	320
27.84.4.2 operator<< . . . . .	320
27.85gdcmm::DictPrinter Class Reference . . . . .	320
27.85.1 Detailed Description . . . . .	322
27.85.2 Constructor & Destructor Documentation . . . . .	322
27.85.2.1 DictPrinter . . . . .	322
27.85.2.2 ~DictPrinter . . . . .	322
27.85.3 Member Function Documentation . . . . .	322
27.85.3.1 Print . . . . .	322
27.85.3.2 PrintDataElement2 . . . . .	322
27.85.3.3 PrintDataSet2 . . . . .	322
27.86gdcmm::Dicts Class Reference . . . . .	322
27.86.1 Detailed Description . . . . .	323
27.86.2 Member Enumeration Documentation . . . . .	323
27.86.2.1 ConstructorType . . . . .	323

27.86.3 Constructor & Destructor Documentation . . . . .	323
27.86.3.1 Dicts . . . . .	323
27.86.3.2 ~Dicts . . . . .	323
27.86.4 Member Function Documentation . . . . .	323
27.86.4.1 GetConstructorString . . . . .	324
27.86.4.2 GetCSAHeaderDict . . . . .	324
27.86.4.3 GetDictEntry . . . . .	324
27.86.4.4 GetDictEntry . . . . .	324
27.86.4.5 GetPrivateDict . . . . .	324
27.86.4.6 GetPrivateDict . . . . .	324
27.86.4.7 GetPublicDict . . . . .	324
27.86.4.8 IsEmpty . . . . .	324
27.86.4.9 LoadDefaults . . . . .	324
27.86.5 Friends And Related Function Documentation . . . . .	324
27.86.5.1 Global . . . . .	324
27.86.5.2 operator<< . . . . .	324
27.87gdcmm::network::DIMSE Class Reference . . . . .	325
27.87.1 Detailed Description . . . . .	325
27.87.2 Member Enumeration Documentation . . . . .	325
27.87.2.1 CommandTypes . . . . .	325
27.88gdcmm::DirectionCosines Class Reference . . . . .	326
27.88.1 Detailed Description . . . . .	327
27.88.2 Constructor & Destructor Documentation . . . . .	327
27.88.2.1 DirectionCosines . . . . .	327
27.88.2.2 DirectionCosines . . . . .	327
27.88.2.3 ~DirectionCosines . . . . .	327
27.88.3 Member Function Documentation . . . . .	327
27.88.3.1 ComputeDistAlongNormal . . . . .	327
27.88.3.2 Cross . . . . .	327
27.88.3.3 CrossDot . . . . .	327
27.88.3.4 Dot . . . . .	327
27.88.3.5 IsValid . . . . .	327
27.88.3.6 Normalize . . . . .	328
27.88.3.7 operator const double * . . . . .	328
27.88.3.8 Print . . . . .	328
27.88.3.9 SetFromString . . . . .	328
27.89gdcmm::Directory Class Reference . . . . .	328

27.89.1 Detailed Description . . . . .	329
27.89.2 Member Typedef Documentation . . . . .	329
27.89.2.1 FilenamesType . . . . .	329
27.89.2.2 FilenameType . . . . .	329
27.89.3 Constructor & Destructor Documentation . . . . .	329
27.89.3.1 Directory . . . . .	329
27.89.3.2 ~Directory . . . . .	329
27.89.4 Member Function Documentation . . . . .	329
27.89.4.1 Explore . . . . .	329
27.89.4.2 GetDirectories . . . . .	330
27.89.4.3 GetFilenames . . . . .	330
27.89.4.4 GetToplevel . . . . .	330
27.89.4.5 Load . . . . .	330
27.89.4.6 Print . . . . .	330
27.89.5 Friends And Related Function Documentation . . . . .	330
27.89.5.1 operator<< . . . . .	330
27.90gdcmm::DirectoryHelper Class Reference . . . . .	331
27.90.1 Detailed Description . . . . .	331
27.90.2 Member Function Documentation . . . . .	331
27.90.2.1 GetCTImageSeriesUIDs . . . . .	331
27.90.2.2 GetFilenamesFromSeriesUIDs . . . . .	331
27.90.2.3 GetFrameOfReference . . . . .	331
27.90.2.4 GetMRImageSeriesUIDs . . . . .	331
27.90.2.5 GetRTStructSeriesUIDs . . . . .	332
27.90.2.6 GetSeriesUIDsBySOPClassUID . . . . .	332
27.90.2.7 GetSOPClassUID . . . . .	332
27.90.2.8 GetStringValueFromTag . . . . .	332
27.90.2.9 LoadImageFromFiles . . . . .	332
27.90.2.10RetrieveSOPInstanceUIDFromIndex . . . . .	332
27.90.2.11RetrieveSOPInstanceUIDFromZPosition . . . . .	332
27.91gdcmm::DummyValueGenerator Class Reference . . . . .	332
27.91.1 Detailed Description . . . . .	332
27.91.2 Member Function Documentation . . . . .	332
27.91.2.1 Generate . . . . .	333
27.92gdcmm::Dumper Class Reference . . . . .	333
27.92.1 Detailed Description . . . . .	334
27.92.2 Constructor & Destructor Documentation . . . . .	334

27.92.2.1 Dumper	334
27.92.2.2 ~Dumper	334
27.93gdcmm::Element< TVR, TVM > Class Template Reference	335
27.93.1 Detailed Description	336
27.93.2 Member Typedef Documentation	337
27.93.2.1 Type	337
27.93.3 Member Function Documentation	337
27.93.3.1 GetAsDataElement	337
27.93.3.2 GetLength	337
27.93.3.3 GetValue	337
27.93.3.4 GetValue	337
27.93.3.5 GetValues	337
27.93.3.6 GetVM	337
27.93.3.7 GetVR	337
27.93.3.8 operator[]	337
27.93.3.9 Print	337
27.93.3.10Read	337
27.93.3.11Set	337
27.93.3.12SetFromDataElement	337
27.93.3.13SetNoSwap	337
27.93.3.14SetValue	337
27.93.3.15Write	337
27.93.4 Member Data Documentation	337
27.93.4.1 Internal	337
27.94gdcmm::Element< TVR, VM::VM1_2 > Class Template Reference	338
27.94.1 Member Typedef Documentation	339
27.94.1.1 Parent	339
27.94.2 Member Function Documentation	339
27.94.2.1 SetLength	339
27.95gdcmm::Element< TVR, VM::VM1_n > Class Template Reference	339
27.95.1 Member Typedef Documentation	340
27.95.1.1 Type	340
27.95.2 Constructor & Destructor Documentation	340
27.95.2.1 Element	340
27.95.2.2 ~Element	340
27.95.2.3 Element	340
27.95.3 Member Function Documentation	340

27.95.3.1 GetAsDataElement . . . . .	340
27.95.3.2 GetLength . . . . .	340
27.95.3.3 GetValue . . . . .	340
27.95.3.4 GetValue . . . . .	340
27.95.3.5 GetVM . . . . .	340
27.95.3.6 GetVR . . . . .	341
27.95.3.7 operator= . . . . .	341
27.95.3.8 operator[] . . . . .	341
27.95.3.9 Print . . . . .	341
27.95.3.10 Read . . . . .	341
27.95.3.11 Set . . . . .	341
27.95.3.12 SetArray . . . . .	341
27.95.3.13 SetFromDataElement . . . . .	341
27.95.3.14 SetLength . . . . .	341
27.95.3.15 SetNoSwap . . . . .	341
27.95.3.16 SetValue . . . . .	341
27.95.3.17 Write . . . . .	341
27.95.3.18 WriteASCII . . . . .	341
27.96gdcm::Element< TVR, VM::VM2_2n > Class Template Reference . . . . .	341
27.96.1 Member Typedef Documentation . . . . .	343
27.96.1.1 Parent . . . . .	343
27.96.2 Member Function Documentation . . . . .	343
27.96.2.1 SetLength . . . . .	343
27.97gdcm::Element< TVR, VM::VM2_n > Class Template Reference . . . . .	343
27.97.1 Member Typedef Documentation . . . . .	344
27.97.1.1 Parent . . . . .	344
27.97.2 Member Function Documentation . . . . .	344
27.97.2.1 SetLength . . . . .	344
27.98gdcm::Element< TVR, VM::VM3_3n > Class Template Reference . . . . .	344
27.98.1 Member Typedef Documentation . . . . .	346
27.98.1.1 Parent . . . . .	346
27.98.2 Member Function Documentation . . . . .	346
27.98.2.1 SetLength . . . . .	346
27.99gdcm::Element< TVR, VM::VM3_n > Class Template Reference . . . . .	346
27.99.1 Member Typedef Documentation . . . . .	347
27.99.1.1 Parent . . . . .	347
27.99.2 Member Function Documentation . . . . .	347

27.99.2.1 SetLength . . . . .	347
27.100dcm::Element< VR::AS, VM::VM5 > Class Template Reference . . . . .	347
27.100.1Member Function Documentation . . . . .	348
27.100.1.1GetLength . . . . .	348
27.100.1.2Print . . . . .	348
27.100.2Member Data Documentation . . . . .	348
27.100.2.1Internal . . . . .	348
27.101dcm::Element< VR::OB, VM::VM1 > Class Template Reference . . . . .	348
27.102dcm::Element< VR::OW, VM::VM1 > Class Template Reference . . . . .	349
27.103dcm::ElementDisableCombinations< TVR, TVM > Class Template Reference . . . . .	351
27.103.1Detailed Description . . . . .	351
27.104dcm::ElementDisableCombinations< VR::OB, VM::VM1_n > Class Template Reference . . . . .	352
27.105dcm::ElementDisableCombinations< VR::OW, VM::VM1_n > Class Template Reference . . . . .	352
27.106dcm::EncapsulatedDocument Class Reference . . . . .	352
27.106.1Detailed Description . . . . .	352
27.106.2Constructor & Destructor Documentation . . . . .	352
27.106.2.1EncapsulatedDocument . . . . .	352
27.107dcm::EncodingImplementation< T > Class Template Reference . . . . .	353
27.107.1Detailed Description . . . . .	353
27.108dcm::EncodingImplementation< VR::VRASCII > Class Template Reference . . . . .	353
27.108.1Member Function Documentation . . . . .	353
27.108.1.1Read . . . . .	353
27.108.1.2ReadComputeLength . . . . .	354
27.108.1.3ReadNoSwap . . . . .	354
27.108.1.4Write . . . . .	354
27.108.1.5Write . . . . .	354
27.108.1.6Write . . . . .	354
27.109dcm::EncodingImplementation< VR::VRBINARY > Class Template Reference . . . . .	354
27.109.1Member Function Documentation . . . . .	354
27.109.1.1Read . . . . .	354
27.109.1.2ReadComputeLength . . . . .	355
27.109.1.3ReadNoSwap . . . . .	355
27.109.1.4Write . . . . .	355
27.110dcm::EndEvent Class Reference . . . . .	355
27.111dcm::EnumeratedValues Class Reference . . . . .	356
27.111.1Detailed Description . . . . .	356
27.111.2Constructor & Destructor Documentation . . . . .	357

27.111.2.1EnumeratedValues . . . . .	357
27.111gdcmm::Event Class Reference . . . . .	357
27.112.1Detailed Description . . . . .	358
27.112.2Constructor & Destructor Documentation . . . . .	358
27.112.2.1Event . . . . .	358
27.112.2.2Event . . . . .	358
27.112.2.3~Event . . . . .	358
27.112.3Member Function Documentation . . . . .	358
27.112.3.1CheckEvent . . . . .	358
27.112.3.2GetEventName . . . . .	358
27.112.3.3MakeObject . . . . .	358
27.112.3.4Print . . . . .	359
27.113gdcmm::Exception Class Reference . . . . .	359
27.113.1Detailed Description . . . . .	360
27.113.2Constructor & Destructor Documentation . . . . .	360
27.113.2.1Exception . . . . .	360
27.113.2.2~Exception . . . . .	360
27.113.3Member Function Documentation . . . . .	360
27.113.3.1GetDescription . . . . .	360
27.113.3.2what . . . . .	360
27.114gdcmm::ExitEvent Class Reference . . . . .	361
27.115gdcmm::ExplicitDataElement Class Reference . . . . .	362
27.115.1Detailed Description . . . . .	363
27.115.2Member Function Documentation . . . . .	363
27.115.2.1GetLength . . . . .	363
27.115.2.2Read . . . . .	363
27.115.2.3ReadPreValue . . . . .	363
27.115.2.4ReadValue . . . . .	363
27.115.2.5ReadWithLength . . . . .	363
27.115.2.6Write . . . . .	363
27.116gdcmm::ExplicitImplicitDataElement Class Reference . . . . .	363
27.116.1Detailed Description . . . . .	365
27.116.2Member Function Documentation . . . . .	365
27.116.2.1GetLength . . . . .	365
27.116.2.2Read . . . . .	365
27.116.2.3ReadPreValue . . . . .	365
27.116.2.4ReadValue . . . . .	365

27.116.2.5	ReadWithLength . . . . .	365
27.117	dcm::Fiducials Class Reference . . . . .	365
27.117.1	Detailed Description . . . . .	365
27.117.2	Constructor & Destructor Documentation . . . . .	365
27.117.2.1	Fiducials . . . . .	365
27.118	dcm::File Class Reference . . . . .	366
27.118.1	Detailed Description . . . . .	367
27.118.2	Constructor & Destructor Documentation . . . . .	367
27.118.2.1	File . . . . .	367
27.118.3	Member Function Documentation . . . . .	367
27.118.3.1	GetDataSet . . . . .	367
27.118.3.2	GetDataSet . . . . .	368
27.118.3.3	GetHeader . . . . .	368
27.118.3.4	GetHeader . . . . .	368
27.118.3.5	Read . . . . .	368
27.118.3.6	SetDataSet . . . . .	368
27.118.3.7	SetHeader . . . . .	368
27.118.3.8	Write . . . . .	368
27.118.4	Friends And Related Function Documentation . . . . .	369
27.118.4.1	operator<< . . . . .	369
27.119	dcm::FileAnonymizer Class Reference . . . . .	369
27.119.1	Detailed Description . . . . .	370
27.119.2	Constructor & Destructor Documentation . . . . .	371
27.119.2.1	FileAnonymizer . . . . .	371
27.119.2.2	~FileAnonymizer . . . . .	371
27.119.3	Member Function Documentation . . . . .	371
27.119.3.1	Empty . . . . .	371
27.119.3.2	Remove . . . . .	371
27.119.3.3	Replace . . . . .	371
27.119.3.4	Replace . . . . .	371
27.119.3.5	SetInputFileName . . . . .	371
27.119.3.6	SetOutputFileName . . . . .	372
27.119.3.7	Write . . . . .	372
27.120	dcm::FileChangeTransferSyntax Class Reference . . . . .	372
27.120.1	Detailed Description . . . . .	373
27.120.2	Constructor & Destructor Documentation . . . . .	374
27.120.2.1	FileChangeTransferSyntax . . . . .	374



27.120.2.2~FileChangeTransferSyntax . . . . .	374
27.120.3 Member Function Documentation . . . . .	374
27.120.3.1Change . . . . .	374
27.120.3.2GetCodec . . . . .	374
27.120.3.3New . . . . .	374
27.120.3.4SetInputFileName . . . . .	374
27.120.3.5SetOutputFileName . . . . .	374
27.120.3.6SetTransferSyntax . . . . .	374
27.121dcm::FileDerivation Class Reference . . . . .	374
27.121.1 Detailed Description . . . . .	375
27.121.2 Constructor & Destructor Documentation . . . . .	375
27.121.2.1FileDerivation . . . . .	375
27.121.2.2~FileDerivation . . . . .	375
27.121.3 Member Function Documentation . . . . .	375
27.121.3.1AddDerivationDescription . . . . .	375
27.121.3.2AddPurposeOfReferenceCodeSequence . . . . .	375
27.121.3.3AddReference . . . . .	376
27.121.3.4AddSourceImageSequence . . . . .	376
27.121.3.5Derive . . . . .	376
27.121.3.6GetFile . . . . .	376
27.121.3.7GetFile . . . . .	376
27.121.3.8SetDerivationCodeSequenceCodeValue . . . . .	376
27.121.3.9SetDerivationDescription . . . . .	376
27.121.3.10SetFile . . . . .	376
27.121.3.11SetPurposeOfReferenceCodeSequenceCodeValue . . . . .	377
27.122dcm::FileExplicitFilter Class Reference . . . . .	377
27.122.1 Detailed Description . . . . .	377
27.122.2 Constructor & Destructor Documentation . . . . .	378
27.122.2.1FileExplicitFilter . . . . .	378
27.122.2.2~FileExplicitFilter . . . . .	378
27.122.3 Member Function Documentation . . . . .	378
27.122.3.1Change . . . . .	378
27.122.3.2ChangeFMI . . . . .	378
27.122.3.3GetFile . . . . .	378
27.122.3.4ProcessDataSet . . . . .	378
27.122.3.5SetChangePrivateTags . . . . .	378
27.122.3.6SetFile . . . . .	378

27.122.3.7SetRecomputeItemLength . . . . .	378
27.122.3.8SetRecomputeSequenceLength . . . . .	379
27.122.3.9SetUseVRUN . . . . .	379
27.123.0dcm::FileMetaInformation Class Reference . . . . .	379
27.123.1Detailed Description . . . . .	381
27.123.2Constructor & Destructor Documentation . . . . .	381
27.123.2.1FileMetaInformation . . . . .	381
27.123.2.2FileMetaInformation . . . . .	381
27.123.3Member Function Documentation . . . . .	381
27.123.3.1AppendImplementationClassUID . . . . .	381
27.123.3.2ComputeDataSetMediaStorageSOPClass . . . . .	381
27.123.3.3ComputeDataSetTransferSyntax . . . . .	381
27.123.3.4Default . . . . .	381
27.123.3.5FillFromDataSet . . . . .	381
27.123.3.6GetDataSetTransferSyntax . . . . .	382
27.123.3.7GetFileMetaInformationVersion . . . . .	382
27.123.3.8GetFullLength . . . . .	382
27.123.3.9GetGDCMImplementationClassUID . . . . .	382
27.123.3.10GetGDCMImplementationVersionName . . . . .	382
27.123.3.11GetGDCMSourceApplicationEntityTitle . . . . .	382
27.123.3.12GetImplementationClassUID . . . . .	382
27.123.3.13GetImplementationVersionName . . . . .	382
27.123.3.14GetMediaStorage . . . . .	382
27.123.3.15GetMediaStorageAsString . . . . .	382
27.123.3.16GetMetaInformationTS . . . . .	382
27.123.3.17GetPreamble . . . . .	382
27.123.3.18GetPreamble . . . . .	382
27.123.3.19GetSourceApplicationEntityTitle . . . . .	382
27.123.3.20Insert . . . . .	382
27.123.3.21IsValid . . . . .	382
27.123.3.22Read . . . . .	382
27.123.3.23ReadCompat . . . . .	383
27.123.3.24ReadCompatInternal . . . . .	383
27.123.3.25Replace . . . . .	383
27.123.3.26SetDataSetTransferSyntax . . . . .	383
27.123.3.27SetImplementationClassUID . . . . .	383
27.123.3.28SetImplementationVersionName . . . . .	383

27.123.3.29	SetPreamble	383
27.123.3.30	SetSourceApplicationEntityTitle	383
27.123.3.31	Write	383
27.123.4	Friends And Related Function Documentation	383
27.123.4.1	operator<<	383
27.123.5	Member Data Documentation	383
27.123.5.1	DataSetMS	383
27.123.5.2	DataSetTS	384
27.123.5.3	MetaInformationTS	384
27.124	dcm::Filename Class Reference	384
27.124.1	Detailed Description	385
27.124.2	Constructor & Destructor Documentation	385
27.124.2.1	Filename	385
27.124.3	Member Function Documentation	385
27.124.3.1	EndWith	385
27.124.3.2	GetExtension	385
27.124.3.3	GetFileName	385
27.124.3.4	GetName	385
27.124.3.5	GetPath	385
27.124.3.6	IsEmpty	385
27.124.3.7	IsIdentical	385
27.124.3.8	Join	385
27.124.3.9	operator const char *	386
27.124.3.10	ToUnixSlashes	386
27.124.3.11	ToWindowsSlashes	386
27.125	dcm::FileNameEvent Class Reference	386
27.125.1	Detailed Description	387
27.125.2	Member Typedef Documentation	388
27.125.2.1	Self	388
27.125.2.2	Superclass	388
27.125.3	Constructor & Destructor Documentation	388
27.125.3.1	FileNameEvent	388
27.125.3.2	~FileNameEvent	388
27.125.3.3	FileNameEvent	388
27.125.4	Member Function Documentation	388
27.125.4.1	CheckEvent	388
27.125.4.2	GetEventName	388

27.125.4.3	GetFileName	388
27.125.4.4	MakeObject	388
27.125.4.5	SetFileName	388
27.126	dcm::FilenameGenerator Class Reference	388
27.126.1	Detailed Description	389
27.126.2	Member Typedef Documentation	389
27.126.2.1	FileNamesType	389
27.126.2.2	FilenameType	389
27.126.2.3	SizeType	389
27.126.3	Constructor & Destructor Documentation	389
27.126.3.1	FilenameGenerator	390
27.126.3.2	~FilenameGenerator	390
27.126.4	Member Function Documentation	390
27.126.4.1	Generate	390
27.126.4.2	GetFilename	390
27.126.4.3	GetFileNames	390
27.126.4.4	GetNumberOfFileNames	390
27.126.4.5	GetPattern	390
27.126.4.6	GetPrefix	390
27.126.4.7	SetNumberOfFileNames	390
27.126.4.8	SetPattern	390
27.126.4.9	SetPrefix	391
27.127	dcm::FileSet Class Reference	391
27.127.1	Detailed Description	391
27.127.2	Member Typedef Documentation	391
27.127.2.1	FilesType	391
27.127.2.2	FileType	391
27.127.3	Constructor & Destructor Documentation	391
27.127.3.1	FileSet	391
27.127.4	Member Function Documentation	391
27.127.4.1	AddFile	392
27.127.4.2	AddFile	392
27.127.4.3	GetFiles	392
27.127.4.4	SetFiles	392
27.127.5	Friends And Related Function Documentation	392
27.127.5.1	operator<<	392
27.128	dcm::FileStreamer Class Reference	392

27.128.1Detailed Description	394
27.128.2Constructor & Destructor Documentation	394
27.128.2.1FileStreamer	394
27.128.2.2~FileStreamer	394
27.128.3Member Function Documentation	394
27.128.3.1AppendToDataElement	394
27.128.3.2AppendToGroupDataElement	394
27.128.3.3CheckDataElement	394
27.128.3.4CheckTemplateFileName	394
27.128.3.5New	394
27.128.3.6ReserveDataElement	395
27.128.3.7ReserveGroupDataElement	395
27.128.3.8SetOutputFileName	395
27.128.3.9SetTemplateFileName	395
27.128.3.10StartDataElement	395
27.128.3.11StartGroupDataElement	395
27.128.3.12StopDataElement	395
27.128.3.13StopGroupDataElement	395
27.129dcm::FileWithName Class Reference	396
27.129.1Detailed Description	397
27.129.2Constructor & Destructor Documentation	397
27.129.2.1FileWithName	397
27.129.3Member Data Documentation	397
27.129.3.1filename	397
27.130dcm::FindPatientRootQuery Class Reference	397
27.130.1Detailed Description	398
27.130.2Constructor & Destructor Documentation	398
27.130.2.1FindPatientRootQuery	398
27.130.3Member Function Documentation	398
27.130.3.1GetAbstractSyntaxUID	398
27.130.3.2GetTagListByLevel	398
27.130.3.3InitializeDataSet	399
27.130.3.4ValidateQuery	399
27.130.4Friends And Related Function Documentation	399
27.130.4.1QueryFactory	399
27.131dcm::FindStudyRootQuery Class Reference	399
27.131.1Detailed Description	401

27.131.2	Constructor & Destructor Documentation	401
27.131.2.1	FindStudyRootQuery	401
27.131.3	Member Function Documentation	401
27.131.3.1	GetAbstractSyntaxUID	401
27.131.3.2	GetTagListByLevel	401
27.131.3.3	InitializeDataSet	401
27.131.3.4	ValidateQuery	401
27.131.4	Friends And Related Function Documentation	401
27.131.4.1	QueryFactory	401
27.132	dcm::Fragment Class Reference	401
27.132.1	Detailed Description	403
27.132.2	Constructor & Destructor Documentation	403
27.132.2.1	Fragment	403
27.132.3	Member Function Documentation	403
27.132.3.1	ComputeLength	403
27.132.3.2	GetLength	403
27.132.3.3	Read	403
27.132.3.4	ReadBacktrack	403
27.132.3.5	ReadPreValue	403
27.132.3.6	ReadValue	403
27.132.3.7	Write	404
27.132.4	Friends And Related Function Documentation	404
27.132.4.1	operator<<	404
27.133	dcm::Global Class Reference	404
27.133.1	Detailed Description	405
27.133.2	Constructor & Destructor Documentation	405
27.133.2.1	Global	405
27.133.2.2	~Global	405
27.133.3	Member Function Documentation	405
27.133.3.1	Append	405
27.133.3.2	GetDefs	405
27.133.3.3	GetDicts	405
27.133.3.4	GetDicts	405
27.133.3.5	GetInstance	405
27.133.3.6	LoadResourcesFiles	406
27.133.3.7	Locate	406
27.133.3.8	Prepend	406

27.133.4	Friends And Related Function Documentation	406
27.133.4.1	operator<<	406
27.134	dcm::GroupDict Class Reference	406
27.134.1	Detailed Description	407
27.134.2	Member Typedef Documentation	407
27.134.2.1	GroupStringVector	407
27.134.3	Constructor & Destructor Documentation	407
27.134.3.1	GroupDict	407
27.134.3.2	~GroupDict	407
27.134.4	Member Function Documentation	407
27.134.4.1	Add	407
27.134.4.2	GetAbbreviation	407
27.134.4.3	GetName	407
27.134.4.4	Insert	408
27.134.4.5	Size	408
27.134.5	Friends And Related Function Documentation	408
27.134.5.1	operator<<	408
27.135	dcm::IconImageFilter Class Reference	408
27.135.1	Detailed Description	408
27.135.2	Constructor & Destructor Documentation	409
27.135.2.1	IconImageFilter	409
27.135.2.2	~IconImageFilter	409
27.135.3	Member Function Documentation	409
27.135.3.1	Extract	409
27.135.3.2	ExtractIconImages	409
27.135.3.3	ExtractVeprolIconImages	409
27.135.3.4	GetFile	409
27.135.3.5	GetFile	409
27.135.3.6	GetIconImage	410
27.135.3.7	GetNumberOfIconImages	410
27.135.3.8	SetFile	410
27.136	dcm::IconImageGenerator Class Reference	410
27.136.1	Detailed Description	411
27.136.2	Constructor & Destructor Documentation	411
27.136.2.1	IconImageGenerator	411
27.136.2.2	~IconImageGenerator	411
27.136.3	Member Function Documentation	411

27.136.3.1AutoPixelMinMax . . . . .	411
27.136.3.2ConvertRGBToPaletteColor . . . . .	411
27.136.3.3Generate . . . . .	411
27.136.3.4GetIconImage . . . . .	412
27.136.3.5GetPixmap . . . . .	412
27.136.3.6GetPixmap . . . . .	412
27.136.3.7SetOutputDimensions . . . . .	412
27.136.3.8SetOutsideValuePixel . . . . .	412
27.136.3.9SetPixelMinMax . . . . .	412
27.136.3.10SetPixmap . . . . .	412
27.137.dcm:ignore_char Struct Reference . . . . .	412
27.137.Constructor & Destructor Documentation . . . . .	413
27.137.1.ignore_char . . . . .	413
27.137.2.Member Data Documentation . . . . .	413
27.137.2.1m_char . . . . .	413
27.138.dcm::Image Class Reference . . . . .	413
27.138.Detailed Description . . . . .	414
27.138.Constructor & Destructor Documentation . . . . .	415
27.138.2.1Image . . . . .	415
27.138.2.2~Image . . . . .	415
27.138.3.Member Function Documentation . . . . .	415
27.138.3.1GetDirectionCosines . . . . .	415
27.138.3.2GetDirectionCosines . . . . .	415
27.138.3.3GetIntercept . . . . .	415
27.138.3.4GetOrigin . . . . .	415
27.138.3.5GetOrigin . . . . .	416
27.138.3.6GetSlope . . . . .	416
27.138.3.7GetSpacing . . . . .	416
27.138.3.8GetSpacing . . . . .	416
27.138.3.9Print . . . . .	416
27.138.3.10SetDirectionCosines . . . . .	416
27.138.3.11SetDirectionCosines . . . . .	416
27.138.3.12SetDirectionCosines . . . . .	416
27.138.3.13SetIntercept . . . . .	416
27.138.3.14SetOrigin . . . . .	416
27.138.3.15SetOrigin . . . . .	416
27.138.3.16SetOrigin . . . . .	416



27.138.3.1SetSlope . . . . .	416
27.138.3.1SetSpacing . . . . .	416
27.138.3.1SetSpacing . . . . .	417
27.139dcm::ImageApplyLookupTable Class Reference . . . . .	417
27.139.1Detailed Description . . . . .	419
27.139.2Constructor & Destructor Documentation . . . . .	419
27.139.2.1ImageApplyLookupTable . . . . .	419
27.139.2.2~ImageApplyLookupTable . . . . .	419
27.139.3Member Function Documentation . . . . .	419
27.139.3.1Apply . . . . .	419
27.140dcm::ImageChangePhotometricInterpretation Class Reference . . . . .	419
27.140.1Detailed Description . . . . .	421
27.140.2Constructor & Destructor Documentation . . . . .	421
27.140.2.1ImageChangePhotometricInterpretation . . . . .	421
27.140.2.2~ImageChangePhotometricInterpretation . . . . .	421
27.140.3Member Function Documentation . . . . .	421
27.140.3.1Change . . . . .	421
27.140.3.2ChangeMonochrome . . . . .	421
27.140.3.3GetPhotometricInterpretation . . . . .	421
27.140.3.4RGB2YBR . . . . .	421
27.140.3.5SetPhotometricInterpretation . . . . .	421
27.140.3.6YBR2RGB . . . . .	422
27.141dcm::ImageChangePlanarConfiguration Class Reference . . . . .	422
27.141.1Detailed Description . . . . .	424
27.141.2Constructor & Destructor Documentation . . . . .	424
27.141.2.1ImageChangePlanarConfiguration . . . . .	424
27.141.2.2~ImageChangePlanarConfiguration . . . . .	424
27.141.3Member Function Documentation . . . . .	424
27.141.3.1Change . . . . .	424
27.141.3.2GetPlanarConfiguration . . . . .	424
27.141.3.3RGBPixelsToRGBPlanes . . . . .	424
27.141.3.4RGBPlanesToRGBPixels . . . . .	424
27.141.3.5SetPlanarConfiguration . . . . .	424
27.142dcm::ImageChangeTransferSyntax Class Reference . . . . .	425
27.142.1Detailed Description . . . . .	427
27.142.2Constructor & Destructor Documentation . . . . .	427
27.142.2.1ImageChangeTransferSyntax . . . . .	427

27.142.2.2~ImageChangeTransferSyntax . . . . .	427
27.142.3Member Function Documentation . . . . .	427
27.142.3.1Change . . . . .	427
27.142.3.2GetTransferSyntax . . . . .	427
27.142.3.3SetCompressIconImage . . . . .	428
27.142.3.4SetForce . . . . .	428
27.142.3.5SetTransferSyntax . . . . .	428
27.142.3.6SetUserCodec . . . . .	428
27.142.3.7TryJPEG2000Codec . . . . .	428
27.142.3.8TryJPEGCodec . . . . .	428
27.142.3.9TryJPEGLSCodec . . . . .	428
27.142.3.10TryRAWCodec . . . . .	428
27.142.3.11TryRLECodec . . . . .	428
27.143dcm::ImageCodec Class Reference . . . . .	429
27.143.1Detailed Description . . . . .	431
27.143.2Member Typedef Documentation . . . . .	431
27.143.2.1LUTPtr . . . . .	431
27.143.3Constructor & Destructor Documentation . . . . .	431
27.143.3.1ImageCodec . . . . .	431
27.143.3.2~ImageCodec . . . . .	431
27.143.4Member Function Documentation . . . . .	431
27.143.4.1AppendFrameEncode . . . . .	431
27.143.4.2AppendRowEncode . . . . .	431
27.143.4.3CanCode . . . . .	432
27.143.4.4CanDecode . . . . .	432
27.143.4.5Clone . . . . .	432
27.143.4.6Decode . . . . .	432
27.143.4.7DecodeByStreams . . . . .	432
27.143.4.8DoByteSwap . . . . .	432
27.143.4.9DoInvertMonochrome . . . . .	432
27.143.4.10DoOverlayCleanup . . . . .	432
27.143.4.11DoPaddedCompositePixelCode . . . . .	432
27.143.4.12DoPlanarConfiguration . . . . .	432
27.143.4.13DoSimpleCopy . . . . .	432
27.143.4.14DoYBR . . . . .	432
27.143.4.15GetDimensions . . . . .	433
27.143.4.16GetHeaderInfo . . . . .	433

27.143.4.10	GetLossyFlag	433
27.143.4.10	GetLUT	433
27.143.4.10	GetNeedByteSwap	433
27.143.4.20	GetNumberOfDimensions	433
27.143.4.20	GetPhotometricInterpretation	433
27.143.4.20	GetPixelFormat	433
27.143.4.20	GetPixelFormat	433
27.143.4.20	GetPlanarConfiguration	433
27.143.4.25	FrameEncoder	433
27.143.4.26	Lossy	433
27.143.4.26	RowEncoder	433
27.143.4.26	Valid	433
27.143.4.29	Dimensions	433
27.143.4.30	Dimensions	434
27.143.4.30	LossyFlag	434
27.143.4.30	LUT	434
27.143.4.30	NeedByteSwap	434
27.143.4.30	NeedOverlayCleanup	434
27.143.4.30	NumberOfDimensions	434
27.143.4.30	PhotometricInterpretation	434
27.143.4.30	PixelFormat	434
27.143.4.30	PlanarConfiguration	434
27.143.4.30	StartEncode	434
27.143.4.40	StopEncode	434
27.143.5	Friends And Related Function Documentation	434
27.143.5.1	FileChangeTransferSyntax	434
27.143.5.2	ImageChangePhotometricInterpretation	435
27.143.6	Member Data Documentation	435
27.143.6.1	Dimensions	435
27.143.6.2	LossyFlag	435
27.143.6.3	LUT	435
27.143.6.4	NeedByteSwap	435
27.143.6.5	NeedOverlayCleanup	435
27.143.6.6	NumberOfDimensions	435
27.143.6.7	PF	435
27.143.6.8	PI	435
27.143.6.9	PlanarConfiguration	435

27.143.6.1	RequestPaddedCompositePixelCode	. . . . .	435
27.143.6.1	RequestPlanarConfiguration	. . . . .	435
27.144	dcm::ImageConverter Class Reference	. . . . .	435
27.144.1	Detailed Description	. . . . .	435
27.144.2	Constructor & Destructor Documentation	. . . . .	436
27.144.2.1	ImageConverter	. . . . .	436
27.144.2.2	~ImageConverter	. . . . .	436
27.144.3	Member Function Documentation	. . . . .	436
27.144.3.1	Convert	. . . . .	436
27.144.3.2	GetOutput	. . . . .	436
27.144.3.3	SetInput	. . . . .	436
27.145	dcm::ImageFragmentSplitter Class Reference	. . . . .	436
27.145.1	Detailed Description	. . . . .	439
27.145.2	Constructor & Destructor Documentation	. . . . .	439
27.145.2.1	ImageFragmentSplitter	. . . . .	439
27.145.2.2	~ImageFragmentSplitter	. . . . .	439
27.145.3	Member Function Documentation	. . . . .	439
27.145.3.1	GetFragmentSizeMax	. . . . .	439
27.145.3.2	SetForce	. . . . .	439
27.145.3.3	SetFragmentSizeMax	. . . . .	439
27.145.3.4	Split	. . . . .	439
27.146	dcm::ImageHelper Class Reference	. . . . .	439
27.146.1	Detailed Description	. . . . .	440
27.146.2	Member Function Documentation	. . . . .	440
27.146.2.1	ComputeMediaStorageFromModality	. . . . .	440
27.146.2.2	ComputeSpacingFromImagePositionPatient	. . . . .	441
27.146.2.3	GetDimensionsValue	. . . . .	441
27.146.2.4	GetDirectionCosinesFromDataSet	. . . . .	441
27.146.2.5	GetDirectionCosinesValue	. . . . .	441
27.146.2.6	GetForcePixelSpacing	. . . . .	441
27.146.2.7	GetForceRescaleInterceptSlope	. . . . .	441
27.146.2.8	GetLUT	. . . . .	441
27.146.2.9	GetOriginValue	. . . . .	441
27.146.2.10	GetPhotometricInterpretationValue	. . . . .	441
27.146.2.11	GetPixelFormatValue	. . . . .	441
27.146.2.12	GetPlanarConfigurationValue	. . . . .	441
27.146.2.13	GetPointerFromElement	. . . . .	441

27.146.2.10	<a href="#">GetRescaleInterceptSlopeValue</a>	441
27.146.2.11	<a href="#">GetSpacingTagFromMediaStorage</a>	442
27.146.2.12	<a href="#">GetSpacingValue</a>	442
27.146.2.13	<a href="#">GetZSpacingTagFromMediaStorage</a>	442
27.146.2.14	<a href="#">SetDimensionsValue</a>	442
27.146.2.15	<a href="#">SetDirectionCosinesValue</a>	442
27.146.2.16	<a href="#">SetForcePixelSpacing</a>	442
27.146.2.17	<a href="#">SetForceRescaleInterceptSlope</a>	442
27.146.2.18	<a href="#">SetOriginValue</a>	442
27.146.2.19	<a href="#">SetRescaleInterceptSlopeValue</a>	442
27.146.2.20	<a href="#">SetSpacingValue</a>	442
27.147	<a href="#">dcm::ImageReader Class Reference</a>	443
27.147.1	<a href="#">Detailed Description</a>	445
27.147.2	<a href="#">Constructor &amp; Destructor Documentation</a>	445
27.147.2.1	<a href="#">ImageReader</a>	445
27.147.2.2	<a href="#">~ImageReader</a>	445
27.147.3	<a href="#">Member Function Documentation</a>	445
27.147.3.1	<a href="#">GetImage</a>	445
27.147.3.2	<a href="#">GetImage</a>	445
27.147.3.3	<a href="#">Read</a>	445
27.147.3.4	<a href="#">ReadACRNEMAImage</a>	446
27.147.3.5	<a href="#">ReadImage</a>	446
27.148	<a href="#">dcm::ImageRegionReader Class Reference</a>	446
27.148.1	<a href="#">Detailed Description</a>	448
27.148.2	<a href="#">Constructor &amp; Destructor Documentation</a>	448
27.148.2.1	<a href="#">ImageRegionReader</a>	448
27.148.2.2	<a href="#">~ImageRegionReader</a>	448
27.148.3	<a href="#">Member Function Documentation</a>	448
27.148.3.1	<a href="#">ComputeBufferLength</a>	448
27.148.3.2	<a href="#">GetRegion</a>	448
27.148.3.3	<a href="#">Read</a>	448
27.148.3.4	<a href="#">ReadInformation</a>	448
27.148.3.5	<a href="#">ReadIntoBuffer</a>	449
27.148.3.6	<a href="#">SetRegion</a>	449
27.149	<a href="#">dcm::ImageToImageFilter Class Reference</a>	449
27.149.1	<a href="#">Detailed Description</a>	450
27.149.2	<a href="#">Constructor &amp; Destructor Documentation</a>	450

27.149.2.1ImageToImageFilter . . . . .	451
27.149.2.2~ImageToImageFilter . . . . .	451
27.149.3Member Function Documentation . . . . .	451
27.149.3.1GetInput . . . . .	451
27.149.3.2GetOutput . . . . .	451
27.150dcm::ImageWriter Class Reference . . . . .	451
27.150.1Detailed Description . . . . .	453
27.150.2Constructor & Destructor Documentation . . . . .	453
27.150.2.1ImageWriter . . . . .	453
27.150.2.2~ImageWriter . . . . .	453
27.150.3Member Function Documentation . . . . .	453
27.150.3.1GetImage . . . . .	453
27.150.3.2GetImage . . . . .	453
27.150.3.3Write . . . . .	453
27.151dcm::network::ImplementationClassUIDSub Class Reference . . . . .	454
27.151.1Detailed Description . . . . .	454
27.151.2Constructor & Destructor Documentation . . . . .	454
27.151.2.1ImplementationClassUIDSub . . . . .	454
27.151.3Member Function Documentation . . . . .	454
27.151.3.1Print . . . . .	454
27.151.3.2Read . . . . .	454
27.151.3.3Size . . . . .	454
27.151.3.4Write . . . . .	454
27.152dcm::network::ImplementationUIDSub Class Reference . . . . .	454
27.152.1Detailed Description . . . . .	455
27.152.2Constructor & Destructor Documentation . . . . .	455
27.152.2.1ImplementationUIDSub . . . . .	455
27.152.3Member Function Documentation . . . . .	455
27.152.3.1Write . . . . .	455
27.153dcm::network::ImplementationVersionNameSub Class Reference . . . . .	455
27.153.1Detailed Description . . . . .	455
27.153.2Constructor & Destructor Documentation . . . . .	455
27.153.2.1ImplementationVersionNameSub . . . . .	455
27.153.3Member Function Documentation . . . . .	455
27.153.3.1Print . . . . .	455
27.153.3.2Read . . . . .	455
27.153.3.3Size . . . . .	455

27.153.3.4Write . . . . .	456
27.154dcm::ImplicitDataElement Class Reference . . . . .	456
27.154.1Detailed Description . . . . .	457
27.154.2Member Function Documentation . . . . .	457
27.154.2.1GetLength . . . . .	457
27.154.2.2Read . . . . .	457
27.154.2.3ReadPreValue . . . . .	457
27.154.2.4ReadValue . . . . .	457
27.154.2.5ReadValueWithLength . . . . .	457
27.154.2.6ReadWithLength . . . . .	457
27.154.2.7Write . . . . .	457
27.155dcm::InitializeEvent Class Reference . . . . .	458
27.156dcm::IOD Class Reference . . . . .	459
27.156.1Detailed Description . . . . .	459
27.156.2Member Typedef Documentation . . . . .	459
27.156.2.1MapIODEntry . . . . .	459
27.156.2.2SizeType . . . . .	459
27.156.3Constructor & Destructor Documentation . . . . .	459
27.156.3.1IOD . . . . .	459
27.156.4Member Function Documentation . . . . .	460
27.156.4.1AddIODEntry . . . . .	460
27.156.4.2Clear . . . . .	460
27.156.4.3GetIODEntry . . . . .	460
27.156.4.4GetNumberOfIODs . . . . .	460
27.156.4.5GetTypeFromTag . . . . .	460
27.156.5Friends And Related Function Documentation . . . . .	460
27.156.5.1operator<< . . . . .	460
27.157dcm::IODEntry Class Reference . . . . .	460
27.157.1Detailed Description . . . . .	461
27.157.2Constructor & Destructor Documentation . . . . .	461
27.157.2.1IODEntry . . . . .	461
27.157.3Member Function Documentation . . . . .	461
27.157.3.1GetIE . . . . .	461
27.157.3.2GetName . . . . .	461
27.157.3.3GetRef . . . . .	461
27.157.3.4GetUsage . . . . .	462
27.157.3.5GetUsageType . . . . .	462

27.157.3.6SetIE	. 462
27.157.3.7SetName	. 462
27.157.3.8SetRef	. 462
27.157.3.9SetUsage	. 462
27.157.4Friends And Related Function Documentation	. 462
27.157.4.1operator<<	. 462
27.158gdcmm::IODs Class Reference	. 462
27.158.1Detailed Description	. 463
27.158.2Member Typedef Documentation	. 463
27.158.2.1IODMapType	. 463
27.158.2.2IODMapTypeConstIterator	. 463
27.158.2.3IODName	. 463
27.158.3Constructor & Destructor Documentation	. 463
27.158.3.1IODs	. 463
27.158.4Member Function Documentation	. 463
27.158.4.1AddIOD	. 463
27.158.4.2Begin	. 463
27.158.4.3Clear	. 463
27.158.4.4End	. 463
27.158.4.5GetIOD	. 463
27.158.5Friends And Related Function Documentation	. 463
27.158.5.1operator<<	. 464
27.159gdcmm::IPPSorter Class Reference	. 464
27.159.1Detailed Description	. 465
27.159.2Constructor & Destructor Documentation	. 465
27.159.2.1IPPSorter	. 465
27.159.3Member Function Documentation	. 465
27.159.3.1GetDirectionCosinesTolerance	. 465
27.159.3.2GetZSpacing	. 466
27.159.3.3GetZSpacingTolerance	. 466
27.159.3.4SetComputeZSpacing	. 466
27.159.3.5SetDirectionCosinesTolerance	. 466
27.159.3.6SetDropDuplicatePositions	. 466
27.159.3.7SetZSpacingTolerance	. 466
27.159.3.8Sort	. 467
27.159.4Member Data Documentation	. 467
27.159.4.1ComputeZSpacing	. 467



27.159.4.2DirCosTolerance . . . . .	467
27.159.4.3DropDuplicatePositions . . . . .	467
27.159.4.4ZSpacing . . . . .	467
27.159.4.5ZTolerance . . . . .	467
27.160dcm::Item Class Reference . . . . .	467
27.160.1Detailed Description . . . . .	468
27.160.2Constructor & Destructor Documentation . . . . .	469
27.160.2.1Item . . . . .	469
27.160.2.2Item . . . . .	469
27.160.3Member Function Documentation . . . . .	469
27.160.3.1Clear . . . . .	469
27.160.3.2FindDataElement . . . . .	469
27.160.3.3GetDataElement . . . . .	469
27.160.3.4GetLength . . . . .	469
27.160.3.5GetNestedDataSet . . . . .	469
27.160.3.6GetNestedDataSet . . . . .	470
27.160.3.7InsertDataElement . . . . .	470
27.160.3.8Read . . . . .	470
27.160.3.9SetNestedDataSet . . . . .	470
27.160.3.10Write . . . . .	470
27.160.4Friends And Related Function Documentation . . . . .	470
27.160.4.1operator<< . . . . .	470
27.161dcm::IterationEvent Class Reference . . . . .	470
27.162dcm::JPEG12Codec Class Reference . . . . .	471
27.162.1Detailed Description . . . . .	473
27.162.2Constructor & Destructor Documentation . . . . .	473
27.162.2.1JPEG12Codec . . . . .	473
27.162.2.2~JPEG12Codec . . . . .	473
27.162.3Member Function Documentation . . . . .	473
27.162.3.1DecodeByStreams . . . . .	473
27.162.3.2EncodeBuffer . . . . .	473
27.162.3.3GetHeaderInfo . . . . .	473
27.162.3.4InternalCode . . . . .	473
27.162.3.5IsStateSuspension . . . . .	473
27.163dcm::JPEG16Codec Class Reference . . . . .	474
27.163.1Detailed Description . . . . .	475
27.163.2Constructor & Destructor Documentation . . . . .	475

27.163.2.1JPEG16Codec . . . . .	475
27.163.2.2~JPEG16Codec . . . . .	475
27.163.3Member Function Documentation . . . . .	475
27.163.3.1DecodeByStreams . . . . .	475
27.163.3.2EncodeBuffer . . . . .	475
27.163.3.3GetHeaderInfo . . . . .	475
27.163.3.4InternalCode . . . . .	475
27.163.3.5IsStateSuspension . . . . .	475
27.164dcm::JPEG2000Codec Class Reference . . . . .	476
27.164.1Detailed Description . . . . .	477
27.164.2Constructor & Destructor Documentation . . . . .	477
27.164.2.1JPEG2000Codec . . . . .	477
27.164.2.2~JPEG2000Codec . . . . .	477
27.164.3Member Function Documentation . . . . .	478
27.164.3.1AppendFrameEncode . . . . .	478
27.164.3.2AppendRowEncode . . . . .	478
27.164.3.3CanCode . . . . .	478
27.164.3.4CanDecode . . . . .	478
27.164.3.5Clone . . . . .	478
27.164.3.6Code . . . . .	478
27.164.3.7Decode . . . . .	478
27.164.3.8DecodeByStreams . . . . .	478
27.164.3.9DecodeExtent . . . . .	478
27.164.3.10GetHeaderInfo . . . . .	479
27.164.3.11GetQuality . . . . .	479
27.164.3.12GetRate . . . . .	479
27.164.3.13FrameEncoder . . . . .	479
27.164.3.14RowEncoder . . . . .	479
27.164.3.15SetNumberOfResolutions . . . . .	479
27.164.3.16SetQuality . . . . .	479
27.164.3.17SetRate . . . . .	479
27.164.3.18SetReversible . . . . .	479
27.164.3.19SetTileSize . . . . .	479
27.164.3.20StartEncode . . . . .	479
27.164.3.21StopEncode . . . . .	479
27.164.4Friends And Related Function Documentation . . . . .	479
27.164.4.1Bitmap . . . . .	479

27.164.4.2	ImageRegionReader	479
27.165	dcm::JPEG8Codec Class Reference	479
27.165.1	Detailed Description	481
27.165.2	Constructor & Destructor Documentation	481
27.165.2.1	JPEG8Codec	481
27.165.2.2	~JPEG8Codec	481
27.165.3	Member Function Documentation	481
27.165.3.1	DecodeByStreams	481
27.165.3.2	EncodeBuffer	481
27.165.3.3	GetHeaderInfo	481
27.165.3.4	InternalCode	481
27.165.3.5	IsStateSuspension	481
27.166	dcm::JPEGCodec Class Reference	482
27.166.1	Detailed Description	483
27.166.2	Constructor & Destructor Documentation	484
27.166.2.1	JPEGCodec	484
27.166.2.2	~JPEGCodec	484
27.166.3	Member Function Documentation	484
27.166.3.1	AppendFrameEncode	484
27.166.3.2	AppendRowEncode	484
27.166.3.3	CanCode	484
27.166.3.4	CanDecode	484
27.166.3.5	Clone	484
27.166.3.6	Code	485
27.166.3.7	ComputeOffsetTable	485
27.166.3.8	Decode	485
27.166.3.9	DecodeByStreams	485
27.166.3.10	DecodeExtent	485
27.166.3.11	EncodeBuffer	485
27.166.3.12	GetHeaderInfo	485
27.166.3.13	GetLossless	485
27.166.3.14	GetQuality	485
27.166.3.15	FrameEncoder	485
27.166.3.16	RowEncoder	485
27.166.3.17	StateSuspension	486
27.166.3.18	Valid	486
27.166.3.19	SetBitSample	486

27.166.3.20	SetLossless	. 486
27.166.3.21	SetPixelFormat	. 486
27.166.3.22	SetQuality	. 486
27.166.3.23	StartEncode	. 486
27.166.3.24	StopEncode	. 486
27.166.4	Friends And Related Function Documentation	. 486
27.166.4.1	ImageRegionReader	. 486
27.166.5	Member Data Documentation	. 486
27.166.5.1	BitSample	. 486
27.166.5.2	Quality	. 486
27.167	dcm::JPEGLSCodec Class Reference	. 486
27.167.1	Detailed Description	. 488
27.167.2	Constructor & Destructor Documentation	. 488
27.167.2.1	JPEGLSCodec	. 488
27.167.2.2	~JPEGLSCodec	. 488
27.167.3	Member Function Documentation	. 488
27.167.3.1	AppendFrameEncode	. 488
27.167.3.2	AppendRowEncode	. 489
27.167.3.3	CanCode	. 489
27.167.3.4	CanDecode	. 489
27.167.3.5	Clone	. 489
27.167.3.6	Code	. 489
27.167.3.7	Decode	. 489
27.167.3.8	Decode	. 489
27.167.3.9	DecodeExtent	. 489
27.167.3.10	GetBufferLength	. 489
27.167.3.11	GetHeaderInfo	. 489
27.167.3.12	GetLossless	. 489
27.167.3.13	FrameEncoder	. 489
27.167.3.14	RowEncoder	. 490
27.167.3.15	SetBufferLength	. 490
27.167.3.16	SetLossless	. 490
27.167.3.17	SetLossyError	. 490
27.167.3.18	StartEncode	. 490
27.167.3.19	StopEncode	. 490
27.167.4	Friends And Related Function Documentation	. 490
27.167.4.1	ImageRegionReader	. 490

27.168	dcm::JSON Class Reference	490
27.168.1	Detailed Description	490
27.168.2	Constructor & Destructor Documentation	491
27.168.2.1	JSON	491
27.168.2.2	~JSON	491
27.168.3	Member Function Documentation	491
27.168.3.1	Code	491
27.168.3.2	Decode	491
27.168.3.3	GetPrettyPrint	491
27.168.3.4	PrettyPrintOff	491
27.168.3.5	PrettyPrintOn	491
27.168.3.6	SetPrettyPrint	491
27.169	dcm::KAKADUCodec Class Reference	491
27.169.1	Detailed Description	493
27.169.2	Constructor & Destructor Documentation	493
27.169.2.1	KAKADUCodec	493
27.169.2.2	~KAKADUCodec	493
27.169.3	Member Function Documentation	493
27.169.3.1	CanCode	493
27.169.3.2	CanDecode	493
27.169.3.3	Clone	493
27.169.3.4	Code	493
27.169.3.5	Decode	493
27.170	dcm::LO Class Reference	493
27.170.1	Detailed Description	495
27.170.2	Member Typedef Documentation	495
27.170.2.1	const_iterator	495
27.170.2.2	const_reference	495
27.170.2.3	const_reverse_iterator	495
27.170.2.4	difference_type	495
27.170.2.5	iterator	495
27.170.2.6	pointer	495
27.170.2.7	reference	495
27.170.2.8	reverse_iterator	495
27.170.2.9	size_type	495
27.170.2.10	Superclass	495
27.170.2.11	Value_type	495

27.170.3	Constructor & Destructor Documentation	495
27.170.3.1	LO	495
27.170.3.2	LO	495
27.170.3.3	LO	495
27.170.3.4	LO	495
27.170.4	Member Function Documentation	495
27.170.4.1	IsValid	495
27.171	gdcm::LookupTable Class Reference	496
27.171.1	Detailed Description	498
27.171.2	Member Enumeration Documentation	498
27.171.2.1	LookupTableType	498
27.171.3	Constructor & Destructor Documentation	498
27.171.3.1	LookupTable	498
27.171.3.2	~LookupTable	498
27.171.3.3	LookupTable	498
27.171.4	Member Function Documentation	498
27.171.4.1	Allocate	498
27.171.4.2	Clear	498
27.171.4.3	Decode	498
27.171.4.4	Decode	499
27.171.4.5	GetBitSample	499
27.171.4.6	GetBufferAsRGBA	499
27.171.4.7	GetLUT	499
27.171.4.8	GetLUTDescriptor	499
27.171.4.9	GetLUTLength	499
27.171.4.10	GetPointer	499
27.171.4.11	InitializeBlueLUT	499
27.171.4.12	Initialized	499
27.171.4.13	InitializeGreenLUT	499
27.171.4.14	InitializeLUT	499
27.171.4.15	InitializeRedLUT	499
27.171.4.16	Print	500
27.171.4.17	SetBlueLUT	500
27.171.4.18	SetGreenLUT	500
27.171.4.19	SetLUT	500
27.171.4.20	SetRedLUT	500
27.171.4.21	WriteBufferAsRGBA	500

27.171.5	Member Data Documentation . . . . .	500
27.171.5.1	BitSample . . . . .	500
27.171.5.2	IncompleteLUT . . . . .	500
27.171.5.3	Internal . . . . .	500
27.172	dcm::Scanner::Itstr Struct Reference . . . . .	500
27.172.1	Member Function Documentation . . . . .	500
27.172.1.1	operator() . . . . .	500
27.173	dcm::Macro Class Reference . . . . .	501
27.173.1	Detailed Description . . . . .	501
27.173.2	Member Typedef Documentation . . . . .	501
27.173.2.1	ArrayIncludeMacrosType . . . . .	501
27.173.2.2	MapModuleEntry . . . . .	501
27.173.3	Constructor & Destructor Documentation . . . . .	501
27.173.3.1	Macro . . . . .	501
27.173.4	Member Function Documentation . . . . .	502
27.173.4.1	AddMacroEntry . . . . .	502
27.173.4.2	Clear . . . . .	502
27.173.4.3	FindMacroEntry . . . . .	502
27.173.4.4	GetMacroEntry . . . . .	502
27.173.4.5	GetName . . . . .	502
27.173.4.6	SetName . . . . .	502
27.173.4.7	Verify . . . . .	502
27.173.5	Friends And Related Function Documentation . . . . .	502
27.173.5.1	operator<< . . . . .	502
27.174	dcm::Macros Class Reference . . . . .	502
27.174.1	Detailed Description . . . . .	503
27.174.2	Member Typedef Documentation . . . . .	503
27.174.2.1	ModuleMapType . . . . .	503
27.174.3	Constructor & Destructor Documentation . . . . .	503
27.174.3.1	Macros . . . . .	503
27.174.4	Member Function Documentation . . . . .	503
27.174.4.1	AddMacro . . . . .	503
27.174.4.2	Clear . . . . .	503
27.174.4.3	GetMacro . . . . .	503
27.174.4.4	IsEmpty . . . . .	503
27.174.5	Friends And Related Function Documentation . . . . .	503
27.174.5.1	operator<< . . . . .	503

27.175	dcm::network::MaximumLengthSub Class Reference . . . . .	503
27.175.1	Detailed Description . . . . .	504
27.175.2	Constructor & Destructor Documentation . . . . .	504
27.175.2.1	MaximumLengthSub . . . . .	504
27.175.3	Member Function Documentation . . . . .	504
27.175.3.1	GetMaximumLength . . . . .	504
27.175.3.2	Print . . . . .	504
27.175.3.3	Read . . . . .	504
27.175.3.4	SetMaximumLength . . . . .	504
27.175.3.5	Size . . . . .	504
27.175.3.6	Write . . . . .	504
27.176	dcm::MD5 Class Reference . . . . .	504
27.176.1	Detailed Description . . . . .	505
27.176.2	Constructor & Destructor Documentation . . . . .	505
27.176.2.1	MD5 . . . . .	505
27.176.2.2	~MD5 . . . . .	505
27.176.3	Member Function Documentation . . . . .	505
27.176.3.1	Compute . . . . .	505
27.176.3.2	ComputeFile . . . . .	505
27.177	dcm::MediaStorage Class Reference . . . . .	505
27.177.1	Detailed Description . . . . .	508
27.177.2	Member Enumeration Documentation . . . . .	508
27.177.2.1	MSType . . . . .	508
27.177.2.2	ObjectType . . . . .	510
27.177.3	Constructor & Destructor Documentation . . . . .	511
27.177.3.1	MediaStorage . . . . .	511
27.177.4	Member Function Documentation . . . . .	511
27.177.4.1	GetModality . . . . .	511
27.177.4.2	GetModalityDimension . . . . .	511
27.177.4.3	GetMSString . . . . .	511
27.177.4.4	GetMSType . . . . .	511
27.177.4.5	GetNumberOfModality . . . . .	511
27.177.4.6	GetNumberOfMSString . . . . .	511
27.177.4.7	GetNumberOfMSType . . . . .	511
27.177.4.8	GetString . . . . .	511
27.177.4.9	GuessFromModality . . . . .	511
27.177.4.10	Image . . . . .	511



27.177.4.1	<del>Is</del> Undefined	512
27.177.4.1	<del>Operator</del> MType	512
27.177.4.1	<del>Set</del> FromDataSet	512
27.177.4.1	<del>Set</del> FromFile	512
27.177.4.1	<del>Set</del> FromHeader	512
27.177.4.1	<del>Set</del> FromModality	512
27.177.4.1	<del>Set</del> FromSourceImageSequence	512
27.177.5	Friends And Related Function Documentation	512
27.177.5.1	<del>operator</del> <<	512
27.178	dcm::MemberCommand< T > Class Template Reference	512
27.178.1	Detailed Description	514
27.178.2	Member Typedef Documentation	514
27.178.2.1	Self	514
27.178.2.2	TConstMemberFunctionPointer	514
27.178.2.3	TMemberFunctionPointer	514
27.178.3	Constructor & Destructor Documentation	515
27.178.3.1	MemberCommand	515
27.178.3.2	~MemberCommand	515
27.178.4	Member Function Documentation	515
27.178.4.1	Execute	515
27.178.4.2	<del>Execute</del>	515
27.178.4.3	New	515
27.178.4.4	SetCallbackFunction	515
27.178.4.5	<del>Set</del> CallbackFunction	515
27.178.5	Member Data Documentation	515
27.178.5.1	m_ConstMemberFunction	516
27.178.5.2	m_MemberFunction	516
27.178.5.3	m_This	516
27.179	dcm::MeshPrimitive Class Reference	516
27.179.1	Detailed Description	518
27.179.2	Member Typedef Documentation	518
27.179.2.1	PrimitivesData	518
27.179.3	Member Enumeration Documentation	518
27.179.3.1	MPType	518
27.179.4	Constructor & Destructor Documentation	519
27.179.4.1	MeshPrimitive	519
27.179.4.2	~MeshPrimitive	519

27.179.5	Member Function Documentation	519
27.179.5.1	AddPrimitiveData	519
27.179.5.2	GetMPType	519
27.179.5.3	GetMPTypeString	519
27.179.5.4	GetNumberOfPrimitivesData	519
27.179.5.5	GetPrimitiveData	519
27.179.5.6	GetPrimitiveData	519
27.179.5.7	GetPrimitiveData	519
27.179.5.8	GetPrimitiveData	519
27.179.5.9	GetPrimitivesData	519
27.179.5.10	GetPrimitivesData	519
27.179.5.11	GetPrimitiveType	519
27.179.5.12	SetPrimitiveData	519
27.179.5.13	SetPrimitiveData	519
27.179.5.14	SetPrimitivesData	519
27.179.5.15	SetPrimitiveType	519
27.179.6	Member Data Documentation	519
27.179.6.1	PrimitiveData	519
27.179.6.2	PrimitiveType	519
27.180	dcm::ModifiedEvent Class Reference	519
27.181	dcm::Module Class Reference	521
27.181.1	Detailed Description	521
27.181.2	Member Typedef Documentation	521
27.181.2.1	ArrayIncludeMacrosType	521
27.181.2.2	MapModuleEntry	521
27.181.3	Constructor & Destructor Documentation	522
27.181.3.1	Module	522
27.181.4	Member Function Documentation	522
27.181.4.1	AddMacro	522
27.181.4.2	AddModuleEntry	522
27.181.4.3	Clear	522
27.181.4.4	FindModuleEntryInMacros	522
27.181.4.5	GetModuleEntryInMacros	522
27.181.4.6	GetName	522
27.181.4.7	SetName	522
27.181.4.8	Verify	522
27.181.5	Friends And Related Function Documentation	522

27.181.5.1operator<<	522
27.182dcm::ModuleEntry Class Reference	522
27.182.1Detailed Description	524
27.182.2Member Typedef Documentation	524
27.182.2.1Description	524
27.182.3Constructor & Destructor Documentation	524
27.182.3.1ModuleEntry	524
27.182.3.2~ModuleEntry	524
27.182.4Member Function Documentation	524
27.182.4.1GetDescription	524
27.182.4.2GetName	524
27.182.4.3GetType	525
27.182.4.4SetDescription	525
27.182.4.5SetName	525
27.182.4.6SetType	525
27.182.5Friends And Related Function Documentation	525
27.182.5.1operator<<	525
27.182.6Member Data Documentation	525
27.182.6.1DataElementType	525
27.182.6.2DescriptionField	525
27.182.6.3Name	525
27.183dcm::Modules Class Reference	525
27.183.1Detailed Description	526
27.183.2Member Typedef Documentation	526
27.183.2.1ModuleMapType	526
27.183.3Constructor & Destructor Documentation	526
27.183.3.1Modules	526
27.183.4Member Function Documentation	526
27.183.4.1AddModule	526
27.183.4.2Clear	526
27.183.4.3GetModule	526
27.183.4.4IsEmpty	526
27.183.5Friends And Related Function Documentation	526
27.183.5.1operator<<	526
27.184dcm::MovePatientRootQuery Class Reference	527
27.184.1Detailed Description	528
27.184.2Constructor & Destructor Documentation	528

27.184.2.1MovePatientRootQuery . . . . .	528
27.184.3Member Function Documentation . . . . .	528
27.184.3.1GetAbstractSyntaxUID . . . . .	528
27.184.3.2GetTagListByLevel . . . . .	528
27.184.3.3InitializeDataSet . . . . .	528
27.184.3.4ValidateQuery . . . . .	528
27.184.4Friends And Related Function Documentation . . . . .	529
27.184.4.1QueryFactory . . . . .	529
27.185dcm::MoveStudyRootQuery Class Reference . . . . .	529
27.185.1Detailed Description . . . . .	530
27.185.2Constructor & Destructor Documentation . . . . .	530
27.185.2.1MoveStudyRootQuery . . . . .	530
27.185.3Member Function Documentation . . . . .	530
27.185.3.1GetAbstractSyntaxUID . . . . .	530
27.185.3.2GetTagListByLevel . . . . .	530
27.185.3.3InitializeDataSet . . . . .	530
27.185.3.4ValidateQuery . . . . .	530
27.185.4Friends And Related Function Documentation . . . . .	531
27.185.4.1QueryFactory . . . . .	531
27.186dcm::NestedModuleEntries Class Reference . . . . .	531
27.186.1Detailed Description . . . . .	532
27.186.2Member Typedef Documentation . . . . .	533
27.186.2.1SizeType . . . . .	533
27.186.3Constructor & Destructor Documentation . . . . .	533
27.186.3.1NestedModuleEntries . . . . .	533
27.186.4Member Function Documentation . . . . .	533
27.186.4.1AddModuleEntry . . . . .	533
27.186.4.2GetModuleEntry . . . . .	533
27.186.4.3GetModuleEntry . . . . .	533
27.186.4.4GetNumberOfModuleEntries . . . . .	533
27.186.5Friends And Related Function Documentation . . . . .	533
27.186.5.1operator<< . . . . .	533
27.187dcm::NoEvent Class Reference . . . . .	533
27.187.1Detailed Description . . . . .	534
27.188dcm::Object Class Reference . . . . .	534
27.188.1Detailed Description . . . . .	536
27.188.2Constructor & Destructor Documentation . . . . .	536

27.188.2.1Object . . . . .	536
27.188.2.2~Object . . . . .	536
27.188.2.3Object . . . . .	536
27.188.3Member Function Documentation . . . . .	536
27.188.3.1operator= . . . . .	536
27.188.3.2Print . . . . .	536
27.188.3.3Register . . . . .	536
27.188.3.4UnRegister . . . . .	536
27.188.4Friends And Related Function Documentation . . . . .	536
27.188.4.1operator<< . . . . .	536
27.188.4.2SmartPointer . . . . .	537
27.189dcm::OpenSSLCryptoFactory Class Reference . . . . .	537
27.189.1Constructor & Destructor Documentation . . . . .	538
27.189.1.1OpenSSLCryptoFactory . . . . .	538
27.189.2Member Function Documentation . . . . .	538
27.189.2.1CreateCMSProvider . . . . .	538
27.189.2.2InitOpenSSL . . . . .	538
27.190dcm::OpenSSLCryptographicMessageSyntax Class Reference . . . . .	538
27.190.1Constructor & Destructor Documentation . . . . .	539
27.190.1.1OpenSSLCryptographicMessageSyntax . . . . .	539
27.190.1.2~OpenSSLCryptographicMessageSyntax . . . . .	539
27.190.2Member Function Documentation . . . . .	539
27.190.2.1Decrypt . . . . .	539
27.190.2.2Encrypt . . . . .	540
27.190.2.3GetCipherType . . . . .	540
27.190.2.4ParseCertificateFile . . . . .	540
27.190.2.5ParseKeyFile . . . . .	540
27.190.2.6SetCipherType . . . . .	540
27.190.2.7SetPassword . . . . .	540
27.191dcm::OpenSSLP7CryptoFactory Class Reference . . . . .	540
27.191.1Constructor & Destructor Documentation . . . . .	541
27.191.1.1OpenSSLP7CryptoFactory . . . . .	541
27.191.2Member Function Documentation . . . . .	541
27.191.2.1CreateCMSProvider . . . . .	542
27.192dcm::OpenSSLP7CryptographicMessageSyntax Class Reference . . . . .	542
27.192.1Detailed Description . . . . .	543
27.192.2Constructor & Destructor Documentation . . . . .	543

27.192.2.1OpenSSLP7CryptographicMessageSyntax . . . . .	543
27.192.2.2~OpenSSLP7CryptographicMessageSyntax . . . . .	543
27.192.3Member Function Documentation . . . . .	543
27.192.3.1Decrypt . . . . .	543
27.192.3.2Encrypt . . . . .	543
27.192.3.3GetCipherType . . . . .	543
27.192.3.4ParseCertificateFile . . . . .	544
27.192.3.5ParseKeyFile . . . . .	544
27.192.3.6SetCipherType . . . . .	544
27.192.3.7SetPassword . . . . .	544
27.193dcm::Orientation Class Reference . . . . .	544
27.193.1Detailed Description . . . . .	545
27.193.2Member Enumeration Documentation . . . . .	545
27.193.2.1OrientationType . . . . .	545
27.193.3Constructor & Destructor Documentation . . . . .	545
27.193.3.1Orientation . . . . .	545
27.193.3.2~Orientation . . . . .	545
27.193.4Member Function Documentation . . . . .	545
27.193.4.1GetLabel . . . . .	545
27.193.4.2GetMajorAxisFromPatientRelativeDirectionCosine . . . . .	546
27.193.4.3GetObliquityThresholdCosineValue . . . . .	546
27.193.4.4GetType . . . . .	546
27.193.4.5Print . . . . .	546
27.193.4.6SetObliquityThresholdCosineValue . . . . .	546
27.193.5Friends And Related Function Documentation . . . . .	546
27.193.5.1operator<< . . . . .	546
27.194dcm::Overlay Class Reference . . . . .	546
27.194.1Detailed Description . . . . .	549
27.194.2Member Enumeration Documentation . . . . .	549
27.194.2.1OverlayType . . . . .	549
27.194.3Constructor & Destructor Documentation . . . . .	549
27.194.3.1Overlay . . . . .	549
27.194.3.2~Overlay . . . . .	549
27.194.3.3Overlay . . . . .	549
27.194.4Member Function Documentation . . . . .	549
27.194.4.1Decompress . . . . .	549
27.194.4.2GetBitPosition . . . . .	549

27.194.4.3	GetBitsAllocated	549
27.194.4.4	GetColumns	549
27.194.4.5	GetDescription	550
27.194.4.6	GetGroup	550
27.194.4.7	GetOrigin	550
27.194.4.8	GetOverlayData	550
27.194.4.9	GetOverlayTypeAsString	550
27.194.4.10	GetOverlayTypeFromString	550
27.194.4.11	GetRows	550
27.194.4.12	GetType	550
27.194.4.13	GetTypeAsEnum	550
27.194.4.14	GetUnpackBuffer	550
27.194.4.15	GetUnpackBufferLength	550
27.194.4.16	GrabOverlayFromPixelData	550
27.194.4.17	IsEmpty	550
27.194.4.18	InPixelData	550
27.194.4.19	InPixelData	551
27.194.4.20	IsZero	551
27.194.4.21	Print	551
27.194.4.22	SetBitPosition	551
27.194.4.23	SetBitsAllocated	551
27.194.4.24	SetColumns	551
27.194.4.25	SetDescription	551
27.194.4.26	SetFrameOrigin	551
27.194.4.27	SetGroup	551
27.194.4.28	SetNumberOfFrames	551
27.194.4.29	SetOrigin	551
27.194.4.30	SetOverlay	552
27.194.4.31	SetRows	552
27.194.4.32	SetType	552
27.194.4.33	Update	552
27.195	dcm::ParseException Class Reference	552
27.195.1	Detailed Description	553
27.195.2	Constructor & Destructor Documentation	553
27.195.2.1	ParseException	553
27.195.2.2	~ParseException	553
27.195.3	Member Function Documentation	553

27.195.3.1	GetLastElement	. . . . .	553
27.195.3.2	operator=	. . . . .	553
27.195.3.3	SetLastElement	. . . . .	553
27.196	dcm::Parser Class Reference	. . . . .	554
27.196.1	Detailed Description	. . . . .	555
27.196.2	Member Typedef Documentation	. . . . .	555
27.196.2.1	EndElementHandler	. . . . .	555
27.196.2.2	StartElementHandler	. . . . .	555
27.196.3	Member Enumeration Documentation	. . . . .	555
27.196.3.1	ErrorType	. . . . .	555
27.196.4	Constructor & Destructor Documentation	. . . . .	555
27.196.4.1	Parser	. . . . .	555
27.196.4.2	~Parser	. . . . .	555
27.196.5	Member Function Documentation	. . . . .	555
27.196.5.1	GetBuffer	. . . . .	555
27.196.5.2	GetCurrentByteIndex	. . . . .	555
27.196.5.3	GetErrorCode	. . . . .	555
27.196.5.4	GetErrorString	. . . . .	555
27.196.5.5	GetUserData	. . . . .	555
27.196.5.6	Parse	. . . . .	556
27.196.5.7	ParseBuffer	. . . . .	556
27.196.5.8	Process	. . . . .	556
27.196.5.9	SetElementHandler	. . . . .	556
27.196.5.10	SetUserData	. . . . .	556
27.197	dcm::Patient Class Reference	. . . . .	556
27.197.1	Detailed Description	. . . . .	556
27.197.2	Constructor & Destructor Documentation	. . . . .	556
27.197.2.1	Patient	. . . . .	556
27.198	dcm::network::PDataTFPDU Class Reference	. . . . .	556
27.198.1	Detailed Description	. . . . .	558
27.198.2	Member Typedef Documentation	. . . . .	558
27.198.2.1	SizeType	. . . . .	558
27.198.3	Constructor & Destructor Documentation	. . . . .	558
27.198.3.1	PDataTFPDU	. . . . .	558
27.198.4	Member Function Documentation	. . . . .	558
27.198.4.1	AddPresentationDataValue	. . . . .	558
27.198.4.2	GetNumberOfPresentationDataValues	. . . . .	558



27.198.4.3	GetPresentationDataValue	558
27.198.4.4	IsLastFragment	558
27.198.4.5	Print	558
27.198.4.6	Read	558
27.198.4.7	ReadInto	558
27.198.4.8	Size	558
27.198.4.9	Write	558
27.199	dcm::PDBElement Class Reference	559
27.199.1	Detailed Description	560
27.199.2	Constructor & Destructor Documentation	560
27.199.2.1	PDBElement	560
27.199.3	Member Function Documentation	560
27.199.3.1	GetName	560
27.199.3.2	GetValue	560
27.199.3.3	operator==	560
27.199.3.4	SetName	560
27.199.3.5	SetValue	560
27.199.4	Friends And Related Function Documentation	560
27.199.4.1	operator<<	560
27.199.5	Member Data Documentation	560
27.199.5.1	NameField	560
27.199.5.2	ValueField	560
27.200	dcm::PDBHeader Class Reference	561
27.200.1	Detailed Description	561
27.200.2	Constructor & Destructor Documentation	562
27.200.2.1	PDBHeader	562
27.200.2.2	~PDBHeader	562
27.200.3	Member Function Documentation	562
27.200.3.1	FindPDBElementByName	562
27.200.3.2	GetPDBEEnd	562
27.200.3.3	GetPDBElementByName	562
27.200.3.4	GetPDBInfoTag	562
27.200.3.5	LoadFromDataElement	562
27.200.3.6	Print	562
27.200.4	Friends And Related Function Documentation	562
27.200.4.1	operator<<	562
27.201	dcm::PDFCodec Class Reference	562

27.201.1	Detailed Description	564
27.201.2	Constructor & Destructor Documentation	564
27.201.2.1	PDFCodec	564
27.201.2.2	PDFCodec	564
27.201.3	Member Function Documentation	564
27.201.3.1	CanCode	564
27.201.3.2	CanDecode	564
27.201.3.3	Decode	564
27.202	dcm::network::PDUFactory Class Reference	564
27.202.1	Detailed Description	565
27.202.2	Member Function Documentation	565
27.202.2.1	ConstructAbortPDU	565
27.202.2.2	ConstructPDU	565
27.202.2.3	ConstructReleasePDU	565
27.202.2.4	CreateCEchoPDU	565
27.202.2.5	CreateCFindPDU	565
27.202.2.6	CreateCMovePDU	565
27.202.2.7	CreateCStoreRQPDU	565
27.202.2.8	CreateCStoreRSPPDU	565
27.202.2.9	DetermineEventByPDU	565
27.202.2.10	GetPDVs	565
27.203	dcm::PersonName Class Reference	566
27.203.1	Detailed Description	566
27.203.2	Member Function Documentation	566
27.203.2.1	GetMaxLength	566
27.203.2.2	GetNumberOfComponents	566
27.203.2.3	Print	566
27.203.2.4	SetBlob	566
27.203.2.5	SetComponents	566
27.203.2.6	SetComponents	566
27.203.3	Member Data Documentation	566
27.203.3.1	Component	567
27.203.3.2	MaxLength	567
27.203.3.3	MaxNumberOfComponents	567
27.203.3.4	Padding	567
27.203.3.5	Separator	567
27.204	dcm::PGXCodec Class Reference	567

27.204.1Detailed Description . . . . .	568
27.204.2Constructor & Destructor Documentation . . . . .	568
27.204.2.1PGXCodec . . . . .	568
27.204.2.2~PGXCodec . . . . .	568
27.204.3Member Function Documentation . . . . .	568
27.204.3.1CanCode . . . . .	568
27.204.3.2CanDecode . . . . .	568
27.204.3.3Clone . . . . .	569
27.204.3.4GetHeaderInfo . . . . .	569
27.204.3.5Read . . . . .	569
27.204.3.6Write . . . . .	569
27.205dcm::PhotometricInterpretation Class Reference . . . . .	569
27.205.1Detailed Description . . . . .	570
27.205.2Member Enumeration Documentation . . . . .	570
27.205.2.1PIType . . . . .	570
27.205.3Constructor & Destructor Documentation . . . . .	571
27.205.3.1PhotometricInterpretation . . . . .	571
27.205.4Member Function Documentation . . . . .	571
27.205.4.1GetPIString . . . . .	571
27.205.4.2GetPIType . . . . .	571
27.205.4.3GetSamplesPerPixel . . . . .	571
27.205.4.4GetString . . . . .	571
27.205.4.5GetType . . . . .	571
27.205.4.6IsLossless . . . . .	571
27.205.4.7IsLossy . . . . .	571
27.205.4.8IsRetired . . . . .	571
27.205.4.9IsSameColorSpace . . . . .	571
27.205.4.10operator PIType . . . . .	571
27.205.5Friends And Related Function Documentation . . . . .	571
27.205.5.1operator<< . . . . .	571
27.206dcm::PixelFormat Class Reference . . . . .	571
27.206.1Detailed Description . . . . .	573
27.206.2Member Enumeration Documentation . . . . .	573
27.206.2.1ScalarType . . . . .	573
27.206.3Constructor & Destructor Documentation . . . . .	574
27.206.3.1PixelFormat . . . . .	574
27.206.3.2PixelFormat . . . . .	574

27.206.4	Member Function Documentation	574
27.206.4.1	GetBitsAllocated	574
27.206.4.2	GetBitsStored	574
27.206.4.3	GetHighBit	574
27.206.4.4	GetMax	574
27.206.4.5	GetMin	575
27.206.4.6	GetPixelRepresentation	575
27.206.4.7	GetPixelSize	575
27.206.4.8	GetSamplesPerPixel	575
27.206.4.9	GetScalarType	575
27.206.4.10	GetScalarTypeAsString	575
27.206.4.11	IsCompatible	575
27.206.4.12	IsValid	575
27.206.4.13	operator ScalarType	575
27.206.4.14	operator !=	575
27.206.4.15	operator !=	575
27.206.4.16	operator ==	576
27.206.4.17	operator ==	576
27.206.4.18	Print	576
27.206.4.19	SetBitsAllocated	576
27.206.4.20	SetBitsStored	576
27.206.4.21	SetHighBit	576
27.206.4.22	SetPixelRepresentation	576
27.206.4.23	SetSamplesPerPixel	576
27.206.4.24	SetScalarType	576
27.206.4.25	Validate	576
27.206.5	Friends And Related Function Documentation	576
27.206.5.1	Bitmap	576
27.206.5.2	operator <<	576
27.207	dcm::Pixmap Class Reference	577
27.207.1	Detailed Description	578
27.207.2	Constructor & Destructor Documentation	578
27.207.2.1	Pixmap	578
27.207.2.2	~Pixmap	578
27.207.3	Member Function Documentation	578
27.207.3.1	AreOverlaysInPixelData	578
27.207.3.2	GetCurve	579

27.207.3.3	GetCurve . . . . .	579
27.207.3.4	GetIconImage . . . . .	579
27.207.3.5	GetIconImage . . . . .	579
27.207.3.6	GetNumberOfCurves . . . . .	579
27.207.3.7	GetNumberOfOverlays . . . . .	579
27.207.3.8	GetOverlay . . . . .	579
27.207.3.9	GetOverlay . . . . .	579
27.207.3.10	Print . . . . .	579
27.207.3.11	RemoveOverlay . . . . .	579
27.207.3.12	SetIconImage . . . . .	579
27.207.3.13	SetNumberOfCurves . . . . .	579
27.207.3.14	SetNumberOfOverlays . . . . .	579
27.207.4	Member Data Documentation . . . . .	579
27.207.4.1	Curves . . . . .	579
27.207.4.2	Icon . . . . .	579
27.207.4.3	Overlays . . . . .	579
27.208	gdcm::PixmapReader Class Reference . . . . .	579
27.208.1	Detailed Description . . . . .	582
27.208.2	Constructor & Destructor Documentation . . . . .	582
27.208.2.1	PixmapReader . . . . .	582
27.208.2.2	~PixmapReader . . . . .	582
27.208.3	Member Function Documentation . . . . .	582
27.208.3.1	GetPixmap . . . . .	582
27.208.3.2	GetPixmap . . . . .	582
27.208.3.3	Read . . . . .	582
27.208.3.4	ReadACRNEMAIImage . . . . .	582
27.208.3.5	ReadImage . . . . .	582
27.208.3.6	ReadImageInternal . . . . .	583
27.208.4	Member Data Documentation . . . . .	583
27.208.4.1	PixelData . . . . .	583
27.209	gdcm::PixmapToPixmapFilter Class Reference . . . . .	583
27.209.1	Detailed Description . . . . .	584
27.209.2	Constructor & Destructor Documentation . . . . .	584
27.209.2.1	PixmapToPixmapFilter . . . . .	584
27.209.2.2	~PixmapToPixmapFilter . . . . .	585
27.209.3	Member Function Documentation . . . . .	585
27.209.3.1	GetInput . . . . .	585

27.209.3.2	GetOutput	. . . . .	585
27.209.3.3	GetOutputAsPixmap	. . . . .	585
27.210	gdcm::PixmapWriter Class Reference	. . . . .	585
27.210.1	Detailed Description	. . . . .	587
27.210.2	Constructor & Destructor Documentation	. . . . .	587
27.210.2.1	PixmapWriter	. . . . .	587
27.210.2.2	~PixmapWriter	. . . . .	587
27.210.3	Member Function Documentation	. . . . .	587
27.210.3.1	IdolconImage	. . . . .	587
27.210.3.2	GetImage	. . . . .	587
27.210.3.3	GetImage	. . . . .	587
27.210.3.4	GetPixmap	. . . . .	587
27.210.3.5	GetPixmap	. . . . .	587
27.210.3.6	PrepareWrite	. . . . .	587
27.210.3.7	SetImage	. . . . .	588
27.210.3.8	SetPixmap	. . . . .	588
27.210.3.9	Write	. . . . .	588
27.210.4	Member Data Documentation	. . . . .	588
27.210.4.1	PixelData	. . . . .	588
27.211	gdcm::PNMCodec Class Reference	. . . . .	588
27.211.1	Detailed Description	. . . . .	590
27.211.2	Constructor & Destructor Documentation	. . . . .	590
27.211.2.1	PNMCodec	. . . . .	590
27.211.2.2	~PNMCodec	. . . . .	590
27.211.3	Member Function Documentation	. . . . .	590
27.211.3.1	CanCode	. . . . .	590
27.211.3.2	CanDecode	. . . . .	590
27.211.3.3	Clone	. . . . .	590
27.211.3.4	GetBufferLength	. . . . .	590
27.211.3.5	GetHeaderInfo	. . . . .	590
27.211.3.6	Read	. . . . .	590
27.211.3.7	SetBufferLength	. . . . .	590
27.211.3.8	Write	. . . . .	591
27.212	gdcm::Preamble Class Reference	. . . . .	591
27.212.1	Detailed Description	. . . . .	591
27.212.2	Constructor & Destructor Documentation	. . . . .	592
27.212.2.1	Preamble	. . . . .	592

27.212.2.2~Preamble . . . . .	592
27.212.2.3Preamble . . . . .	592
27.212.3Member Function Documentation . . . . .	592
27.212.3.1Clear . . . . .	592
27.212.3.2Create . . . . .	592
27.212.3.3GetInternal . . . . .	592
27.212.3.4GetLength . . . . .	592
27.212.3.5IsEmpty . . . . .	592
27.212.3.6IsValid . . . . .	592
27.212.3.7operator= . . . . .	592
27.212.3.8Print . . . . .	592
27.212.3.9Read . . . . .	592
27.212.3.10Remove . . . . .	592
27.212.3.11Valid . . . . .	592
27.212.3.12Write . . . . .	592
27.212.4Friends And Related Function Documentation . . . . .	592
27.212.4.1operator<< . . . . .	592
27.213gdcmm::PresentationContext Class Reference . . . . .	592
27.213.1Detailed Description . . . . .	593
27.213.2Member Typedef Documentation . . . . .	593
27.213.2.1SizeType . . . . .	593
27.213.2.2TransferSyntaxArrayType . . . . .	593
27.213.3Constructor & Destructor Documentation . . . . .	593
27.213.3.1PresentationContext . . . . .	593
27.213.3.2PresentationContext . . . . .	593
27.213.4Member Function Documentation . . . . .	593
27.213.4.1AddTransferSyntax . . . . .	593
27.213.4.2GetAbstractSyntax . . . . .	593
27.213.4.3GetNumberOfTransferSyntaxes . . . . .	593
27.213.4.4GetPresentationContextID . . . . .	594
27.213.4.5GetTransferSyntax . . . . .	594
27.213.4.6operator== . . . . .	594
27.213.4.7Print . . . . .	594
27.213.4.8SetAbstractSyntax . . . . .	594
27.213.4.9SetPresentationContextID . . . . .	594
27.214gdcmm::network::PresentationContextAC Class Reference . . . . .	594
27.214.1Detailed Description . . . . .	594

27.214.2	Constructor & Destructor Documentation	594
27.214.2.1	PresentationContextAC	594
27.214.3	Member Function Documentation	594
27.214.3.1	GetPresentationContextID	595
27.214.3.2	GetReason	595
27.214.3.3	GetTransferSyntax	595
27.214.3.4	Print	595
27.214.3.5	Read	595
27.214.3.6	SetPresentationContextID	595
27.214.3.7	SetReason	595
27.214.3.8	SetTransferSyntax	595
27.214.3.9	Size	595
27.214.3.10	Write	595
27.215	dcm::PresentationContextGenerator Class Reference	595
27.215.1	Detailed Description	596
27.215.2	Member Typedef Documentation	596
27.215.2.1	PresentationContextArrayType	596
27.215.2.2	SizeType	596
27.215.3	Constructor & Destructor Documentation	596
27.215.3.1	PresentationContextGenerator	596
27.215.4	Member Function Documentation	596
27.215.4.1	AddPresentationContext	596
27.215.4.2	GenerateFromFilenames	596
27.215.4.3	GenerateFromUID	597
27.215.4.4	GetDefaultTransferSyntax	597
27.215.4.5	GetPresentationContexts	597
27.215.4.6	SetDefaultTransferSyntax	597
27.215.4.7	SetMergeModeToAbstractSyntax	597
27.215.4.8	SetMergeModeToTransferSyntax	597
27.216	dcm::network::PresentationContextRQ Class Reference	597
27.216.1	Detailed Description	598
27.216.2	Member Typedef Documentation	598
27.216.2.1	SizeType	598
27.216.3	Constructor & Destructor Documentation	598
27.216.3.1	PresentationContextRQ	598
27.216.3.2	PresentationContextRQ	598
27.216.3.3	PresentationContextRQ	598



27.216.4	Member Function Documentation	598
27.216.4.1	AddTransferSyntax	598
27.216.4.2	GetAbstractSyntax	598
27.216.4.3	GetAbstractSyntax	598
27.216.4.4	GetNumberOfTransferSyntaxes	598
27.216.4.5	GetPresentationContextID	598
27.216.4.6	GetTransferSyntax	599
27.216.4.7	GetTransferSyntax	599
27.216.4.8	GetTransferSyntaxes	599
27.216.4.9	operator==	599
27.216.4.10	Print	599
27.216.4.11	Read	599
27.216.4.12	SetAbstractSyntax	599
27.216.4.13	SetPresentationContextID	599
27.216.4.14	Size	599
27.216.4.15	Write	599
27.217	gdcm::network::PresentationDataValue Class Reference	599
27.217.1	Detailed Description	600
27.217.2	Constructor & Destructor Documentation	600
27.217.2.1	PresentationDataValue	600
27.217.3	Member Function Documentation	600
27.217.3.1	ConcatenatePDVBlobs	600
27.217.3.2	ConcatenatePDVBlobsAsExplicit	600
27.217.3.3	GetBlob	600
27.217.3.4	GetIsCommand	600
27.217.3.5	GetIsLastFragment	600
27.217.3.6	GetMessageHeader	600
27.217.3.7	GetPresentationContextID	600
27.217.3.8	Print	600
27.217.3.9	Read	600
27.217.3.10	ReadInto	600
27.217.3.11	SetBlob	600
27.217.3.12	SetCommand	600
27.217.3.13	SetDataSet	600
27.217.3.14	SetLastFragment	601
27.217.3.15	SetMessageHeader	601
27.217.3.16	SetPresentationContextID	601

27.217.3.1Size	601
27.217.3.1Write	601
27.218dcm::Printer Class Reference	601
27.218.1Detailed Description	603
27.218.2Member Enumeration Documentation	603
27.218.2.1PrintStyles	603
27.218.3Constructor & Destructor Documentation	603
27.218.3.1Printer	603
27.218.3.2~Printer	603
27.218.4Member Function Documentation	603
27.218.4.1GetPrintStyle	603
27.218.4.2Print	603
27.218.4.3PrintDataElement	603
27.218.4.4PrintDataSet	603
27.218.4.5PrintSQ	604
27.218.4.6SetColor	604
27.218.4.7SetFile	604
27.218.4.8SetStyle	604
27.218.5Member Data Documentation	604
27.218.5.1F	604
27.218.5.2MaxPrintLength	604
27.218.5.3PrintStyle	604
27.219dcm::PrivateDict Class Reference	604
27.219.1Detailed Description	605
27.219.2Constructor & Destructor Documentation	605
27.219.2.1PrivateDict	605
27.219.2.2~PrivateDict	605
27.219.3Member Function Documentation	605
27.219.3.1AddDictEntry	605
27.219.3.2FindDictEntry	605
27.219.3.3GetDictEntry	605
27.219.3.4IsEmpty	605
27.219.3.5LoadDefault	605
27.219.3.6PrintXML	605
27.219.3.7RemoveDictEntry	605
27.219.4Friends And Related Function Documentation	605
27.219.4.1Dicts	605

27.219.4.2operator<<	605
27.220dcm::PrivateTag Class Reference	606
27.220.1Detailed Description	607
27.220.2Constructor & Destructor Documentation	607
27.220.2.1PrivateTag	607
27.220.2.2PrivateTag	607
27.220.3Member Function Documentation	607
27.220.3.1GetAsDataElement	607
27.220.3.2GetOwner	607
27.220.3.3operator<	607
27.220.3.4ReadFromCommaSeparatedString	607
27.220.3.5SetOwner	607
27.220.4Friends And Related Function Documentation	607
27.220.4.1operator<<	607
27.221dcm::ProgressEvent Class Reference	608
27.221.1Detailed Description	609
27.221.2Member Typedef Documentation	609
27.221.2.1Self	609
27.221.2.2Superclass	609
27.221.3Constructor & Destructor Documentation	609
27.221.3.1ProgressEvent	609
27.221.3.2~ProgressEvent	609
27.221.3.3ProgressEvent	609
27.221.4Member Function Documentation	609
27.221.4.1CheckEvent	609
27.221.4.2GetEventName	609
27.221.4.3GetProgress	609
27.221.4.4MakeObject	609
27.221.4.5SetProgress	610
27.222dcm::PVRGCodec Class Reference	610
27.222.1Detailed Description	611
27.222.2Constructor & Destructor Documentation	611
27.222.2.1PVRGCodec	611
27.222.2.2~PVRGCodec	611
27.222.3Member Function Documentation	611
27.222.3.1CanCode	611
27.222.3.2CanDecode	611

27.222.3.3Clone	. . . . .	612
27.222.3.4Code	. . . . .	612
27.222.3.5Decode	. . . . .	612
27.222.3.6SetLossyFlag	. . . . .	612
27.223dcm::PythonFilter Class Reference	. . . . .	612
27.223.1Detailed Description	. . . . .	612
27.223.2Constructor & Destructor Documentation	. . . . .	613
27.223.2.1PythonFilter	. . . . .	613
27.223.2.2~PythonFilter	. . . . .	613
27.223.3Member Function Documentation	. . . . .	613
27.223.3.1GetFile	. . . . .	613
27.223.3.2GetFile	. . . . .	613
27.223.3.3SetDicts	. . . . .	613
27.223.3.4SetFile	. . . . .	613
27.223.3.5ToPyObject	. . . . .	613
27.223.3.6UseDictAlways	. . . . .	613
27.224dcm::QueryBase Class Reference	. . . . .	613
27.224.1Detailed Description	. . . . .	614
27.224.2Constructor & Destructor Documentation	. . . . .	614
27.224.2.1~QueryBase	. . . . .	614
27.224.3Member Function Documentation	. . . . .	614
27.224.3.1GetAllRequiredTags	. . . . .	614
27.224.3.2GetAllTags	. . . . .	614
27.224.3.3GetHierarchicalSearchTags	. . . . .	614
27.224.3.4GetName	. . . . .	615
27.224.3.5GetOptionalTags	. . . . .	615
27.224.3.6GetQueryLevel	. . . . .	615
27.224.3.7GetRequiredTags	. . . . .	615
27.224.3.8GetUniqueTags	. . . . .	615
27.225dcm::QueryFactory Class Reference	. . . . .	615
27.225.1Detailed Description	. . . . .	615
27.225.2Member Function Documentation	. . . . .	616
27.225.2.1GetCharacterFromCurrentLocale	. . . . .	616
27.225.2.2ListCharSets	. . . . .	616
27.225.2.3ProduceCharacterSetDataElement	. . . . .	616
27.225.2.4ProduceQuery	. . . . .	616
27.226dcm::QueryImage Class Reference	. . . . .	616

27.226.1Detailed Description . . . . .	617
27.226.2Member Function Documentation . . . . .	617
27.226.2.1GetHierarchicalSearchTags . . . . .	618
27.226.2.2GetName . . . . .	618
27.226.2.3GetOptionalTags . . . . .	618
27.226.2.4GetQueryLevel . . . . .	618
27.226.2.5GetRequiredTags . . . . .	618
27.226.2.6GetUniqueTags . . . . .	618
27.227dcm::QueryPatient Class Reference . . . . .	618
27.227.1Detailed Description . . . . .	619
27.227.2Member Function Documentation . . . . .	619
27.227.2.1GetHierarchicalSearchTags . . . . .	620
27.227.2.2GetName . . . . .	620
27.227.2.3GetOptionalTags . . . . .	620
27.227.2.4GetQueryLevel . . . . .	620
27.227.2.5GetRequiredTags . . . . .	620
27.227.2.6GetUniqueTags . . . . .	620
27.228dcm::QuerySeries Class Reference . . . . .	620
27.228.1Detailed Description . . . . .	621
27.228.2Member Function Documentation . . . . .	621
27.228.2.1GetHierarchicalSearchTags . . . . .	622
27.228.2.2GetName . . . . .	622
27.228.2.3GetOptionalTags . . . . .	622
27.228.2.4GetQueryLevel . . . . .	622
27.228.2.5GetRequiredTags . . . . .	622
27.228.2.6GetUniqueTags . . . . .	622
27.229dcm::QueryStudy Class Reference . . . . .	622
27.229.1Detailed Description . . . . .	623
27.229.2Member Function Documentation . . . . .	623
27.229.2.1GetHierarchicalSearchTags . . . . .	624
27.229.2.2GetName . . . . .	624
27.229.2.3GetOptionalTags . . . . .	624
27.229.2.4GetQueryLevel . . . . .	624
27.229.2.5GetRequiredTags . . . . .	624
27.229.2.6GetUniqueTags . . . . .	624
27.230dcm::RAWCodec Class Reference . . . . .	624
27.230.1Detailed Description . . . . .	626

27.230.2	Constructor & Destructor Documentation	626
27.230.2.1	RAWCodec	626
27.230.2.2	~RAWCodec	626
27.230.3	Member Function Documentation	626
27.230.3.1	CanCode	626
27.230.3.2	CanDecode	626
27.230.3.3	Clone	626
27.230.3.4	Code	626
27.230.3.5	Decode	626
27.230.3.6	DecodeByStreams	627
27.230.3.7	DecodeBytes	627
27.230.3.8	GetHeaderInfo	627
27.230.4	gdcm::Reader Class Reference	627
27.231.1	Detailed Description	629
27.231.2	Constructor & Destructor Documentation	630
27.231.2.1	Reader	630
27.231.2.2	~Reader	630
27.231.3	Member Function Documentation	630
27.231.3.1	CanRead	630
27.231.3.2	GetFile	630
27.231.3.3	GetFile	630
27.231.3.4	GetStreamCurrentPosition	630
27.231.3.5	GetStreamPtr	630
27.231.3.6	Read	630
27.231.3.7	ReadDataSet	631
27.231.3.8	ReadMetaInformation	631
27.231.3.9	ReadPreamble	631
27.231.3.10	ReadSelectedPrivateTags	631
27.231.3.11	ReadSelectedTags	631
27.231.3.12	ReadUpToTag	631
27.231.3.13	SetFile	631
27.231.3.14	SetFileName	631
27.231.3.15	SetStream	632
27.231.4	Friends And Related Function Documentation	632
27.231.4.1	StreamImageReader	632
27.231.5	Member Data Documentation	632
27.231.5.1	F	632

27.232	dcm::Region Class Reference . . . . .	632
27.232.1	Detailed Description . . . . .	633
27.232.2	Constructor & Destructor Documentation . . . . .	633
27.232.2.1	Region . . . . .	633
27.232.2.2	~Region . . . . .	633
27.232.3	Member Function Documentation . . . . .	633
27.232.3.1	Area . . . . .	633
27.232.3.2	Clone . . . . .	633
27.232.3.3	ComputeBoundingBox . . . . .	633
27.232.3.4	Empty . . . . .	633
27.232.3.5	IsValid . . . . .	634
27.232.3.6	Print . . . . .	634
27.233	dcm::Rescaler Class Reference . . . . .	634
27.233.1	Detailed Description . . . . .	635
27.233.2	Constructor & Destructor Documentation . . . . .	635
27.233.2.1	Rescaler . . . . .	635
27.233.2.2	~Rescaler . . . . .	635
27.233.3	Member Function Documentation . . . . .	635
27.233.3.1	ComputeInterceptSlopePixelType . . . . .	635
27.233.3.2	ComputePixelTypeFromMinMax . . . . .	636
27.233.3.3	GetIntercept . . . . .	636
27.233.3.4	GetSlope . . . . .	636
27.233.3.5	InverseRescale . . . . .	636
27.233.3.6	InverseRescaleFunctionIntoBestFit . . . . .	636
27.233.3.7	Rescale . . . . .	636
27.233.3.8	RescaleFunctionIntoBestFit . . . . .	636
27.233.3.9	SetIntercept . . . . .	636
27.233.3.10	SetMinMaxForPixelType . . . . .	636
27.233.3.11	SetPixelFormat . . . . .	636
27.233.3.12	SetSlope . . . . .	636
27.233.3.13	SetTargetPixelType . . . . .	636
27.233.3.14	SetUseTargetPixelType . . . . .	636
27.234	dcm::RLECodec Class Reference . . . . .	637
27.234.1	Detailed Description . . . . .	638
27.234.2	Constructor & Destructor Documentation . . . . .	638
27.234.2.1	RLECodec . . . . .	638
27.234.2.2	~RLECodec . . . . .	638

27.234.3	Member Function Documentation	638
27.234.3.1	AppendFrameEncode	639
27.234.3.2	AppendRowEncode	639
27.234.3.3	CanCode	639
27.234.3.4	CanDecode	639
27.234.3.5	Clone	639
27.234.3.6	Code	639
27.234.3.7	Decode	639
27.234.3.8	DecodeByStreams	639
27.234.3.9	DecodeExtent	639
27.234.3.10	GetBufferLength	639
27.234.3.11	GetHeaderInfo	639
27.234.3.12	FrameEncoder	640
27.234.3.13	RowEncoder	640
27.234.3.14	SetBufferLength	640
27.234.3.15	SetLength	640
27.234.3.16	StartEncode	640
27.234.3.17	StopEncode	640
27.234.4	Friends And Related Function Documentation	640
27.234.4.1	ImageRegionReader	640
27.235	dcm::network::RoleSelectionSub Class Reference	640
27.235.1	Detailed Description	640
27.235.2	Constructor & Destructor Documentation	641
27.235.2.1	RoleSelectionSub	641
27.235.3	Member Function Documentation	641
27.235.3.1	Print	641
27.235.3.2	Read	641
27.235.3.3	SetTuple	641
27.235.3.4	Size	641
27.235.3.5	Write	641
27.236	dcm::SerieHelper::Rule Struct Reference	641
27.236.1	Member Data Documentation	642
27.236.1.1	elem	642
27.236.1.2	group	642
27.236.1.3	op	642
27.236.1.4	value	642
27.237	dcm::Scanner Class Reference	642



27.237.1Detailed Description . . . . .	644
27.237.2Member Typedef Documentation . . . . .	645
27.237.2.1ConstIterator . . . . .	645
27.237.2.2MappingType . . . . .	645
27.237.2.3TagToValue . . . . .	645
27.237.2.4TagToValueValueType . . . . .	645
27.237.2.5ValuesType . . . . .	645
27.237.3Constructor & Destructor Documentation . . . . .	645
27.237.3.1Scanner . . . . .	645
27.237.3.2~Scanner . . . . .	645
27.237.4Member Function Documentation . . . . .	645
27.237.4.1AddPrivateTag . . . . .	645
27.237.4.2AddSkipTag . . . . .	645
27.237.4.3AddTag . . . . .	645
27.237.4.4Begin . . . . .	646
27.237.4.5ClearSkipTags . . . . .	646
27.237.4.6ClearTags . . . . .	646
27.237.4.7End . . . . .	646
27.237.4.8GetAllFileNamesFromTagToValue . . . . .	646
27.237.4.9GetFilenameFromTagToValue . . . . .	646
27.237.4.10GetFileNames . . . . .	646
27.237.4.11GetKeys . . . . .	646
27.237.4.12GetMapping . . . . .	646
27.237.4.13GetMappingFromTagToValue . . . . .	646
27.237.4.14GetMappings . . . . .	646
27.237.4.15GetOrderedValues . . . . .	646
27.237.4.16GetValue . . . . .	647
27.237.4.17GetValues . . . . .	647
27.237.4.18GetValues . . . . .	647
27.237.4.19IsKey . . . . .	647
27.237.4.20New . . . . .	647
27.237.4.21Print . . . . .	647
27.237.4.22ProcessPublicTag . . . . .	647
27.237.4.23Scan . . . . .	647
27.237.5Friends And Related Function Documentation . . . . .	648
27.237.5.1operator<< . . . . .	648
27.237.6dcm::Segment Class Reference . . . . .	648

27.238.1Detailed Description . . . . .	650
27.238.2Member Typedef Documentation . . . . .	650
27.238.2.1SurfaceVector . . . . .	650
27.238.3Member Enumeration Documentation . . . . .	650
27.238.3.1ALGOType . . . . .	650
27.238.4Constructor & Destructor Documentation . . . . .	650
27.238.4.1Segment . . . . .	650
27.238.4.2~Segment . . . . .	650
27.238.5Member Function Documentation . . . . .	650
27.238.5.1AddSurface . . . . .	650
27.238.5.2GetALGOType . . . . .	650
27.238.5.3GetALGOTypeString . . . . .	650
27.238.5.4GetAnatomicRegion . . . . .	650
27.238.5.5GetAnatomicRegion . . . . .	650
27.238.5.6GetPropertyCategory . . . . .	650
27.238.5.7GetPropertyCategory . . . . .	650
27.238.5.8GetPropertyType . . . . .	650
27.238.5.9GetPropertyType . . . . .	650
27.238.5.10GetSegmentAlgorithmName . . . . .	651
27.238.5.11GetSegmentAlgorithmType . . . . .	651
27.238.5.12GetSegmentDescription . . . . .	651
27.238.5.13GetSegmentLabel . . . . .	651
27.238.5.14GetSegmentNumber . . . . .	651
27.238.5.15GetSurface . . . . .	651
27.238.5.16GetSurfaceCount . . . . .	651
27.238.5.17GetSurfaces . . . . .	651
27.238.5.18GetSurfaces . . . . .	651
27.238.5.19SetAnatomicRegion . . . . .	651
27.238.5.20SetPropertyCategory . . . . .	651
27.238.5.21SetPropertyType . . . . .	651
27.238.5.22SetSegmentAlgorithmName . . . . .	651
27.238.5.23SetSegmentAlgorithmType . . . . .	651
27.238.5.24SetSegmentAlgorithmType . . . . .	651
27.238.5.25SetSegmentDescription . . . . .	651
27.238.5.26SetSegmentLabel . . . . .	651
27.238.5.27SetSegmentNumber . . . . .	651
27.238.5.28SetSurfaceCount . . . . .	651

27.238.6	Member Data Documentation	651
27.238.6.1	AnatomicRegion	651
27.238.6.2	PropertyCategory	651
27.238.6.3	PropertyType	651
27.238.6.4	SegmentAlgorithmName	651
27.238.6.5	SegmentAlgorithmType	651
27.238.6.6	SegmentDescription	651
27.238.6.7	SegmentLabel	652
27.238.6.8	SegmentNumber	652
27.238.6.9	SurfaceCount	652
27.238.6.10	Surfaces	652
27.239	dcm::SegmentedPaletteColorLookupTable Class Reference	652
27.239.1	Detailed Description	653
27.239.2	Constructor & Destructor Documentation	653
27.239.2.1	SegmentedPaletteColorLookupTable	653
27.239.2.2	~SegmentedPaletteColorLookupTable	653
27.239.3	Member Function Documentation	653
27.239.3.1	Print	653
27.239.3.2	SetLUT	654
27.240	dcm::SegmentReader Class Reference	654
27.240.1	Detailed Description	656
27.240.2	Member Typedef Documentation	656
27.240.2.1	SegmentMap	656
27.240.2.2	SegmentVector	656
27.240.3	Constructor & Destructor Documentation	656
27.240.3.1	SegmentReader	656
27.240.3.2	~SegmentReader	656
27.240.4	Member Function Documentation	656
27.240.4.1	GetSegments	656
27.240.4.2	GetSegments	656
27.240.4.3	Read	656
27.240.4.4	ReadSegment	656
27.240.4.5	ReadSegments	656
27.240.5	Member Data Documentation	656
27.240.5.1	Segments	656
27.241	dcm::SegmentWriter Class Reference	656
27.241.1	Detailed Description	658

27.241.2	Member Typedef Documentation	658
27.241.2.1	SegmentVector	658
27.241.3	Constructor & Destructor Documentation	658
27.241.3.1	SegmentWriter	658
27.241.3.2	~SegmentWriter	658
27.241.4	Member Function Documentation	658
27.241.4.1	AddSegment	658
27.241.4.2	GetNumberOfSegments	658
27.241.4.3	GetSegment	658
27.241.4.4	GetSegments	658
27.241.4.5	GetSegments	658
27.241.4.6	PrepareWrite	658
27.241.4.7	SetNumberOfSegments	658
27.241.4.8	SetSegments	658
27.241.4.9	Write	658
27.241.5	Member Data Documentation	659
27.241.5.1	Segments	659
27.242	gdcmm::SequenceOfFragments Class Reference	659
27.242.1	Detailed Description	661
27.242.2	Member Typedef Documentation	661
27.242.2.1	ConstIterator	661
27.242.2.2	FragmentVector	661
27.242.2.3	Iterator	661
27.242.2.4	SizeType	661
27.242.3	Constructor & Destructor Documentation	661
27.242.3.1	SequenceOfFragments	661
27.242.4	Member Function Documentation	662
27.242.4.1	AddFragment	662
27.242.4.2	Begin	662
27.242.4.3	Begin	662
27.242.4.4	Clear	662
27.242.4.5	ComputeByteLength	662
27.242.4.6	ComputeLength	662
27.242.4.7	End	662
27.242.4.8	End	662
27.242.4.9	GetBuffer	662
27.242.4.10	GetFragBuffer	662

27.242.4.1	<del>Get</del> Fragment	662
27.242.4.1	<del>Get</del> Length	662
27.242.4.1	<del>Get</del> NumberOfFragments	662
27.242.4.1	<del>Get</del> Table	663
27.242.4.1	<del>Get</del> Table	663
27.242.4.1	<del>New</del>	663
27.242.4.1	<del>operator</del> ==	663
27.242.4.1	<del>Print</del>	663
27.242.4.1	<del>Read</del>	663
27.242.4.2	<del>Read</del> PreValue	663
27.242.4.2	<del>Read</del> Value	663
27.242.4.2	<del>Set</del> Length	663
27.242.4.2	<del>Write</del>	663
27.242.4.2	<del>Write</del> Buffer	663
27.243	gdcmm::SequenceOfItems Class Reference	664
27.243.1	Detailed Description	666
27.243.2	Member Typedef Documentation	667
27.243.2.1	ConstIterator	667
27.243.2.2	ItemVector	667
27.243.2.3	Iterator	667
27.243.2.4	SizeType	667
27.243.3	Constructor & Destructor Documentation	667
27.243.3.1	SequenceOfItems	667
27.243.4	Member Function Documentation	667
27.243.4.1	AddItem	667
27.243.4.2	Begin	667
27.243.4.3	Begin	667
27.243.4.4	Clear	667
27.243.4.5	ComputeLength	667
27.243.4.6	End	667
27.243.4.7	End	668
27.243.4.8	FindDataElement	668
27.243.4.9	GetItem	668
27.243.4.10	GetItem	668
27.243.4.10	GetLength	668
27.243.4.10	GetNumberOfItems	668
27.243.4.11	UndefinedLength	668

27.243.4.1New	. 668
27.243.4.15operator=	. 668
27.243.4.16operator==	. 668
27.243.4.17Print	. 669
27.243.4.18Read	. 669
27.243.4.19RemoveItemByIndex	. 669
27.243.4.20SetLength	. 669
27.243.4.21SetLengthToUndefined	. 669
27.243.4.22SetNumberOfItems	. 669
27.243.4.23Write	. 669
27.243.5Member Data Documentation	. 669
27.243.5.1Items	. 669
27.243.5.2SequenceLengthField	. 670
27.244gdcmm::SerieHelper Class Reference	. 670
27.244.1Detailed Description	. 671
27.244.2Member Typedef Documentation	. 671
27.244.2.1SerieRestrictions	. 671
27.244.2.2SingleSerieUIDFileSetmap	. 671
27.244.3Constructor & Destructor Documentation	. 671
27.244.3.1SerieHelper	. 671
27.244.3.2~SerieHelper	. 672
27.244.4Member Function Documentation	. 672
27.244.4.1AddFile	. 672
27.244.4.2AddFileName	. 672
27.244.4.3AddRestriction	. 672
27.244.4.4AddRestriction	. 672
27.244.4.5AddRestriction	. 672
27.244.4.6Clear	. 672
27.244.4.7CreateDefaultUniqueSeriesIdentifier	. 672
27.244.4.8CreateUniqueSeriesIdentifier	. 672
27.244.4.9FileNameOrdering	. 672
27.244.4.10GetFirstSingleSerieUIDFileSet	. 672
27.244.4.11GetNextSingleSerieUIDFileSet	. 672
27.244.4.12ImagePositionPatientOrdering	. 672
27.244.4.13OrderFileList	. 672
27.244.4.14SetDirectory	. 672
27.244.4.15SetLoadMode	. 672

27.244.4.1	<del>Set</del> UseSeriesDetails . . . . .	672
27.244.4.1	<del>User</del> Ordering . . . . .	672
27.244.5	Member Data Documentation . . . . .	672
27.244.5.1	ItFileSetHt . . . . .	672
27.244.5.2	SingleSerieUIDFileSetHT . . . . .	672
27.245	dcm::Series Class Reference . . . . .	672
27.245.1	Detailed Description . . . . .	673
27.245.2	Constructor & Destructor Documentation . . . . .	673
27.245.2.1	Series . . . . .	673
27.246	dcm::network::ServiceClassApplicationInformation Class Reference . . . . .	673
27.246.1	Detailed Description . . . . .	673
27.246.2	Constructor & Destructor Documentation . . . . .	673
27.246.2.1	ServiceClassApplicationInformation . . . . .	673
27.246.3	Member Function Documentation . . . . .	673
27.246.3.1	Print . . . . .	673
27.246.3.2	Read . . . . .	673
27.246.3.3	SetTuple . . . . .	674
27.246.3.4	Size . . . . .	674
27.246.3.5	Write . . . . .	674
27.247	dcm::ServiceClassUser Class Reference . . . . .	674
27.247.1	Detailed Description . . . . .	676
27.247.2	Constructor & Destructor Documentation . . . . .	676
27.247.2.1	ServiceClassUser . . . . .	676
27.247.2.2	~ServiceClassUser . . . . .	676
27.247.3	Member Function Documentation . . . . .	676
27.247.3.1	GetAETitle . . . . .	676
27.247.3.2	GetCalledAETitle . . . . .	677
27.247.3.3	GetTimeout . . . . .	677
27.247.3.4	InitializeConnection . . . . .	677
27.247.3.5	IsPresentationContextAccepted . . . . .	677
27.247.3.6	New . . . . .	677
27.247.3.7	SendEcho . . . . .	677
27.247.3.8	SendFind . . . . .	677
27.247.3.9	SendMove . . . . .	677
27.247.3.10	SendMove . . . . .	677
27.247.3.11	SendMove . . . . .	677
27.247.3.12	SendStore . . . . .	677

27.247.3.1SendStore . . . . .	678
27.247.3.1SendStore . . . . .	678
27.247.3.1SetAETitle . . . . .	678
27.247.3.1SetCalledAETitle . . . . .	678
27.247.3.1SetHostname . . . . .	678
27.247.3.1SetPort . . . . .	678
27.247.3.1SetPortSCP . . . . .	678
27.247.3.2SetPresentationContexts . . . . .	678
27.247.3.2SetTimeout . . . . .	679
27.247.3.2StartAssociation . . . . .	679
27.247.3.2StopAssociation . . . . .	679
27.248dcm::SHA1 Class Reference . . . . .	679
27.248.1Detailed Description . . . . .	680
27.248.2Constructor & Destructor Documentation . . . . .	680
27.248.2.1SHA1 . . . . .	680
27.248.2.2~SHA1 . . . . .	680
27.248.3Member Function Documentation . . . . .	680
27.248.3.1Compute . . . . .	680
27.248.3.2ComputeFile . . . . .	680
27.249dcm::SimpleMemberCommand< T > Class Template Reference . . . . .	680
27.249.1Detailed Description . . . . .	682
27.249.2Member Typedef Documentation . . . . .	682
27.249.2.1Self . . . . .	682
27.249.2.2TMemberFunctionPointer . . . . .	682
27.249.3Constructor & Destructor Documentation . . . . .	683
27.249.3.1SimpleMemberCommand . . . . .	683
27.249.3.2~SimpleMemberCommand . . . . .	683
27.249.4Member Function Documentation . . . . .	683
27.249.4.1Execute . . . . .	683
27.249.4.2Execute . . . . .	683
27.249.4.3New . . . . .	683
27.249.4.4SetCallbackFunction . . . . .	683
27.249.5Member Data Documentation . . . . .	683
27.249.5.1m_MemberFunction . . . . .	683
27.249.5.2m_This . . . . .	684
27.250dcm::SimpleSubjectWatcher Class Reference . . . . .	684
27.250.1Detailed Description . . . . .	684



27.250.2	Constructor & Destructor Documentation	684
27.250.2.1	SimpleSubjectWatcher	684
27.250.2.2	~SimpleSubjectWatcher	684
27.250.3	Member Function Documentation	684
27.250.3.1	EndFilter	685
27.250.3.2	ShowAbort	685
27.250.3.3	ShowAnonymization	685
27.250.3.4	ShowData	685
27.250.3.5	ShowDataSet	685
27.250.3.6	ShowFileName	685
27.250.3.7	ShowIteration	685
27.250.3.8	ShowProgress	685
27.250.3.9	StartFilter	685
27.250.3.10	TestAbortOff	685
27.250.3.11	TestAbortOn	685
27.251	gdcmm::SmartPointer< ObjectType > Class Template Reference	685
27.251.1	Detailed Description	687
27.251.2	Constructor & Destructor Documentation	687
27.251.2.1	SmartPointer	687
27.251.2.2	SmartPointer	687
27.251.2.3	SmartPointer	687
27.251.2.4	SmartPointer	687
27.251.2.5	~SmartPointer	687
27.251.3	Member Function Documentation	687
27.251.3.1	GetPointer	687
27.251.3.2	operator ObjectType *	688
27.251.3.3	operator*	688
27.251.3.4	operator->	688
27.251.3.5	operator=	688
27.251.3.6	operator=	688
27.251.3.7	operator=	688
27.252	gdcmm::network::SOPClassExtendedNegociationSub Class Reference	688
27.252.1	Detailed Description	689
27.252.2	Constructor & Destructor Documentation	689
27.252.2.1	ISOPClassExtendedNegociationSub	689
27.252.3	Member Function Documentation	689
27.252.3.1	Print	689

27.252.3.2	Read	689
27.252.3.3	SetTuple	689
27.252.3.4	Size	689
27.252.3.5	Write	689
27.253	dcm::SOPClassUIDToIOD Class Reference	689
27.253.1	Detailed Description	690
27.253.2	Member Typedef Documentation	690
27.253.2.1	const	690
27.253.3	Member Function Documentation	690
27.253.3.1	GetIOD	690
27.253.3.2	GetIODFromSOPClassUID	690
27.253.3.3	GetNumberOfSOPClassToIOD	690
27.253.3.4	GetSOPClassUIDFromIOD	690
27.253.3.5	GetSOPClassUIDToIOD	690
27.253.3.6	GetSOPClassUIDToIODs	690
27.254	dcm::Sorter Class Reference	690
27.254.1	Detailed Description	692
27.254.2	Member Typedef Documentation	692
27.254.2.1	SelectionMap	692
27.254.2.2	SortFunction	692
27.254.3	Constructor & Destructor Documentation	693
27.254.3.1	Sorter	693
27.254.3.2	~Sorter	693
27.254.4	Member Function Documentation	693
27.254.4.1	AddSelect	693
27.254.4.2	GetFileNames	693
27.254.4.3	Print	693
27.254.4.4	SetSortFunction	693
27.254.4.5	Sort	693
27.254.4.6	StableSort	693
27.254.5	Friends And Related Function Documentation	694
27.254.5.1	operator<<	694
27.254.6	Member Data Documentation	694
27.254.6.1	FileNames	694
27.254.6.2	Selection	694
27.254.6.3	SortFunc	694
27.255	dcm::Spacing Class Reference	694

27.255.1Detailed Description . . . . .	694
27.255.2Member Enumeration Documentation . . . . .	695
27.255.2.1SpacingType . . . . .	695
27.255.3Constructor & Destructor Documentation . . . . .	695
27.255.3.1Spacing . . . . .	695
27.255.3.2~Spacing . . . . .	695
27.255.4Member Function Documentation . . . . .	695
27.255.4.1ComputePixelAspectRatioFromPixelSpacing . . . . .	695
27.256gdcmm::Spectroscopy Class Reference . . . . .	696
27.256.1Detailed Description . . . . .	696
27.256.2Constructor & Destructor Documentation . . . . .	696
27.256.2.1Spectroscopy . . . . .	696
27.257gdcmm::SplitMosaicFilter Class Reference . . . . .	696
27.257.1Detailed Description . . . . .	697
27.257.2Constructor & Destructor Documentation . . . . .	697
27.257.2.1SplitMosaicFilter . . . . .	697
27.257.2.2~SplitMosaicFilter . . . . .	697
27.257.3Member Function Documentation . . . . .	697
27.257.3.1ComputeMOSAICDimensions . . . . .	697
27.257.3.2GetFile . . . . .	697
27.257.3.3GetFile . . . . .	697
27.257.3.4GetImage . . . . .	697
27.257.3.5GetImage . . . . .	697
27.257.3.6SetFile . . . . .	697
27.257.3.7SetImage . . . . .	697
27.257.3.8Split . . . . .	697
27.258gdcmm::StartEvent Class Reference . . . . .	697
27.259gdcmm::static_assert_test< x > Struct Template Reference . . . . .	699
27.260gdcmm::STATIC_ASSERTION_FAILURE< x > Struct Template Reference . . . . .	699
27.261gdcmm::STATIC_ASSERTION_FAILURE< true > Struct Template Reference . . . . .	699
27.261.1Member Enumeration Documentation . . . . .	699
27.261.1.1anonymous enum . . . . .	699
27.262gdcmm::StreamImageReader Class Reference . . . . .	699
27.262.1Detailed Description . . . . .	700
27.262.2Constructor & Destructor Documentation . . . . .	700
27.262.2.1StreamImageReader . . . . .	700
27.262.2.2~StreamImageReader . . . . .	700

27.262.3	Member Function Documentation	700
27.262.3.1	CanReadImage	700
27.262.3.2	DefinePixelExtent	700
27.262.3.3	DefineProperBufferLength	701
27.262.3.4	GetDimensionsValueForResolution	701
27.262.3.5	GetFile	701
27.262.3.6	Read	701
27.262.3.7	ReadImageInformation	701
27.262.3.8	SetFileName	702
27.262.3.9	SetStream	702
27.263	dcm::StreamImageWriter Class Reference	702
27.263.1	Detailed Description	704
27.263.2	Constructor & Destructor Documentation	704
27.263.2.1	StreamImageWriter	704
27.263.2.2	~StreamImageWriter	704
27.263.3	Member Function Documentation	704
27.263.3.1	CanWriteFile	705
27.263.3.2	DefinePixelExtent	705
27.263.3.3	DefineProperBufferLength	705
27.263.3.4	SetFile	705
27.263.3.5	SetFileName	705
27.263.3.6	SetStream	705
27.263.3.7	Write	706
27.263.3.8	WriteImageInformation	706
27.263.3.9	WriteImageSubregionRAW	706
27.263.3.10	WriteRawHeader	706
27.263.4	Member Data Documentation	706
27.263.4.1	mElementOffsets	706
27.263.4.2	mElementOffsets1	707
27.263.4.3	mospFile	707
27.263.4.4	mWriter	707
27.263.4.5	mXMax	707
27.263.4.6	mXMin	707
27.263.4.7	mYMax	707
27.263.4.8	mYMin	707
27.263.4.9	mZMax	707
27.263.4.10	mZMin	707

27.264	dcm::String< TDelimiter, TMaxLength, TPadChar > Class Template Reference . . . . .	707
27.264.1	Detailed Description . . . . .	709
27.264.2	Member Typedef Documentation . . . . .	709
27.264.2.1	const_iterator . . . . .	709
27.264.2.2	const_reference . . . . .	709
27.264.2.3	const_reverse_iterator . . . . .	709
27.264.2.4	difference_type . . . . .	709
27.264.2.5	iterator . . . . .	709
27.264.2.6	pointer . . . . .	709
27.264.2.7	reference . . . . .	709
27.264.2.8	reverse_iterator . . . . .	709
27.264.2.9	size_type . . . . .	709
27.264.2.10	value_type . . . . .	709
27.264.3	Constructor & Destructor Documentation . . . . .	709
27.264.3.1	String . . . . .	710
27.264.3.2	String . . . . .	710
27.264.3.3	String . . . . .	710
27.264.3.4	String . . . . .	710
27.264.4	Member Function Documentation . . . . .	710
27.264.4.1	IsValid . . . . .	710
27.264.4.2	operator const char * . . . . .	710
27.264.4.3	Trim . . . . .	710
27.264.4.4	Trim . . . . .	710
27.264.4.5	Truncate . . . . .	710
27.265	dcm::StringFilter Class Reference . . . . .	710
27.265.1	Detailed Description . . . . .	711
27.265.2	Constructor & Destructor Documentation . . . . .	711
27.265.2.1	StringFilter . . . . .	711
27.265.2.2	~StringFilter . . . . .	711
27.265.3	Member Function Documentation . . . . .	711
27.265.3.1	ExecuteQuery . . . . .	711
27.265.3.2	ExecuteQuery . . . . .	712
27.265.3.3	FromString . . . . .	712
27.265.3.4	FromString . . . . .	712
27.265.3.5	GetFile . . . . .	712
27.265.3.6	GetFile . . . . .	712
27.265.3.7	SetDicts . . . . .	712

27.265.3.8SetFile . . . . .	712
27.265.3.9ToString . . . . .	712
27.265.3.10ToString . . . . .	712
27.265.3.11ToStringPair . . . . .	712
27.265.3.12ToStringPair . . . . .	712
27.265.3.13ToStringPair . . . . .	713
27.265.3.14UseDictAlways . . . . .	713
27.266dcm::Study Class Reference . . . . .	713
27.266.1Detailed Description . . . . .	713
27.266.2Constructor & Destructor Documentation . . . . .	713
27.266.2.1Study . . . . .	713
27.267dcm::Subject Class Reference . . . . .	713
27.267.1Detailed Description . . . . .	715
27.267.2Constructor & Destructor Documentation . . . . .	715
27.267.2.1Subject . . . . .	715
27.267.2.2~Subject . . . . .	715
27.267.3Member Function Documentation . . . . .	715
27.267.3.1AddObserver . . . . .	715
27.267.3.2AddObserver . . . . .	715
27.267.3.3GetCommand . . . . .	715
27.267.3.4HasObserver . . . . .	715
27.267.3.5InvokeEvent . . . . .	715
27.267.3.6InvokeEvent . . . . .	715
27.267.3.7RemoveAllObservers . . . . .	716
27.267.3.8RemoveObserver . . . . .	716
27.268dcm::Surface Class Reference . . . . .	716
27.268.1Detailed Description . . . . .	718
27.268.2Member Enumeration Documentation . . . . .	719
27.268.2.1STATES . . . . .	719
27.268.2.2VIEWType . . . . .	719
27.268.3Constructor & Destructor Documentation . . . . .	719
27.268.3.1Surface . . . . .	719
27.268.3.2~Surface . . . . .	719
27.268.4Member Function Documentation . . . . .	719
27.268.4.1GetAlgorithmFamily . . . . .	719
27.268.4.2GetAlgorithmFamily . . . . .	719
27.268.4.3GetAlgorithmName . . . . .	719

27.268.4.4GetAlgorithmVersion . . . . .	719
27.268.4.5GetAxisOfRotation . . . . .	719
27.268.4.6GetCenterOfRotation . . . . .	720
27.268.4.7GetFiniteVolume . . . . .	720
27.268.4.8GetManifold . . . . .	720
27.268.4.9GetMaximumPointDistance . . . . .	720
27.268.4.10GetMeanPointDistance . . . . .	720
27.268.4.10GetMeshPrimitive . . . . .	720
27.268.4.10GetMeshPrimitive . . . . .	720
27.268.4.10GetNumberOfSurfacePoints . . . . .	720
27.268.4.10GetNumberOfVectors . . . . .	720
27.268.4.10GetPointCoordinatesData . . . . .	720
27.268.4.10GetPointCoordinatesData . . . . .	720
27.268.4.10GetPointPositionAccuracy . . . . .	720
27.268.4.10GetPointsBoundingBoxCoordinates . . . . .	720
27.268.4.10GetProcessingAlgorithm . . . . .	720
27.268.4.20GetProcessingAlgorithm . . . . .	720
27.268.4.20GetRecommendedDisplayCIELabValue . . . . .	720
27.268.4.20GetRecommendedDisplayCIELabValue . . . . .	720
27.268.4.20GetRecommendedDisplayGrayscaleValue . . . . .	720
27.268.4.20GetRecommendedPresentationOpacity . . . . .	720
27.268.4.25GetRecommendedPresentationType . . . . .	721
27.268.4.26GetSTATES . . . . .	721
27.268.4.27GetSTATESString . . . . .	721
27.268.4.28GetSurfaceComments . . . . .	721
27.268.4.29GetSurfaceNumber . . . . .	721
27.268.4.30GetSurfaceProcessing . . . . .	721
27.268.4.30GetSurfaceProcessingDescription . . . . .	721
27.268.4.30GetSurfaceProcessingRatio . . . . .	721
27.268.4.30GetVectorAccuracy . . . . .	721
27.268.4.30GetVectorCoordinateData . . . . .	721
27.268.4.30GetVectorCoordinateData . . . . .	721
27.268.4.30GetVectorDimensionality . . . . .	721
27.268.4.30GetVIEWType . . . . .	721
27.268.4.38GetVIEWTypeString . . . . .	721
27.268.4.39GetAlgorithmFamily . . . . .	721
27.268.4.40GetAlgorithmName . . . . .	721

27.268.4.49	SetAlgorithmVersion	721
27.268.4.49	SetAxisOfRotation	721
27.268.4.49	SetCenterOfRotation	721
27.268.4.49	SetFiniteVolume	721
27.268.4.49	SetManifold	721
27.268.4.49	SetMaximumPointDistance	721
27.268.4.49	SetMeanPointDistance	721
27.268.4.49	SetMeshPrimitive	721
27.268.4.49	SetNumberOfSurfacePoints	721
27.268.4.50	SetNumberOfVectors	721
27.268.4.53	SetPointCoordinatesData	722
27.268.4.53	SetPointPositionAccuracy	722
27.268.4.53	SetPointsBoundingBoxCoordinates	722
27.268.4.53	SetProcessingAlgorithm	722
27.268.4.55	SetRecommendedDisplayCIELabValue	722
27.268.4.55	SetRecommendedDisplayCIELabValue	722
27.268.4.55	SetRecommendedDisplayCIELabValue	722
27.268.4.58	SetRecommendedDisplayGrayscaleValue	722
27.268.4.59	SetRecommendedPresentationOpacity	722
27.268.4.60	SetRecommendedPresentationType	722
27.268.4.63	SetSurfaceComments	722
27.268.4.63	SetSurfaceNumber	722
27.268.4.63	SetSurfaceProcessing	722
27.268.4.63	SetSurfaceProcessingDescription	722
27.268.4.65	SetSurfaceProcessingRatio	722
27.268.4.66	SetVectorAccuracy	722
27.268.4.67	SetVectorCoordinateData	722
27.268.4.68	SetVectorDimensionality	722
27.269	gdcm::SurfaceHelper Class Reference	722
27.269.1	Detailed Description	723
27.269.2	Member Typedef Documentation	723
27.269.2.1	ColorArray	723
27.269.3	Member Function Documentation	723
27.269.3.1	RecommendedDisplayCIELabToRGB	723
27.269.3.2	RecommendedDisplayCIELabToRGB	723
27.269.3.3	RGBToRecommendedDisplayCIELab	724
27.269.3.4	RGBToRecommendedDisplayGrayscale	724



27.270	<a href="#">gdcm::SurfaceReader Class Reference</a>	. 725
27.270.1	<a href="#">Detailed Description</a>	. 726
27.270.2	<a href="#">Constructor &amp; Destructor Documentation</a>	. 727
27.270.2.1	<a href="#">SurfaceReader</a>	. 727
27.270.2.2	<a href="#">~SurfaceReader</a>	. 727
27.270.3	<a href="#">Member Function Documentation</a>	. 727
27.270.3.1	<a href="#">GetNumberOfSurfaces</a>	. 727
27.270.3.2	<a href="#">Read</a>	. 727
27.270.3.3	<a href="#">ReadPointMacro</a>	. 727
27.270.3.4	<a href="#">ReadSurface</a>	. 727
27.270.3.5	<a href="#">ReadSurfaces</a>	. 727
27.271	<a href="#">gdcm::SurfaceWriter Class Reference</a>	. 727
27.271.1	<a href="#">Detailed Description</a>	. 729
27.271.2	<a href="#">Constructor &amp; Destructor Documentation</a>	. 729
27.271.2.1	<a href="#">SurfaceWriter</a>	. 729
27.271.2.2	<a href="#">~SurfaceWriter</a>	. 729
27.271.3	<a href="#">Member Function Documentation</a>	. 729
27.271.3.1	<a href="#">ComputeNumberOfSurfaces</a>	. 729
27.271.3.2	<a href="#">GetNumberOfSurfaces</a>	. 729
27.271.3.3	<a href="#">PrepareWrite</a>	. 729
27.271.3.4	<a href="#">PrepareWritePointMacro</a>	. 729
27.271.3.5	<a href="#">SetNumberOfSurfaces</a>	. 729
27.271.3.6	<a href="#">Write</a>	. 729
27.271.4	<a href="#">Member Data Documentation</a>	. 729
27.271.4.1	<a href="#">NumberOfSurfaces</a>	. 729
27.272	<a href="#">gdcm::SwapCode Class Reference</a>	. 729
27.272.1	<a href="#">Detailed Description</a>	. 730
27.272.2	<a href="#">Member Enumeration Documentation</a>	. 730
27.272.2.1	<a href="#">SwapCodeType</a>	. 730
27.272.3	<a href="#">Constructor &amp; Destructor Documentation</a>	. 731
27.272.3.1	<a href="#">SwapCode</a>	. 731
27.272.4	<a href="#">Member Function Documentation</a>	. 731
27.272.4.1	<a href="#">GetIndex</a>	. 731
27.272.4.2	<a href="#">GetSwapCodeString</a>	. 731
27.272.4.3	<a href="#">operator SwapCode::SwapCodeType</a>	. 731
27.272.5	<a href="#">Friends And Related Function Documentation</a>	. 731
27.272.5.1	<a href="#">operator&lt;&lt;</a>	. 731

27.273	dcm::SwapperDoOp Class Reference . . . . .	731
27.273.1	Member Function Documentation . . . . .	731
27.273.1.1	Swap . . . . .	731
27.273.1.2	SwapArray . . . . .	731
27.274	dcm::SwapperNoOp Class Reference . . . . .	732
27.274.1	Detailed Description . . . . .	732
27.274.2	Member Function Documentation . . . . .	732
27.274.2.1	Swap . . . . .	732
27.274.2.2	SwapArray . . . . .	732
27.275	dcm::System Class Reference . . . . .	732
27.275.1	Detailed Description . . . . .	733
27.275.2	Member Function Documentation . . . . .	733
27.275.2.1	DeleteDirectory . . . . .	733
27.275.2.2	EncodeBytes . . . . .	733
27.275.2.3	FileExists . . . . .	734
27.275.2.4	FileIsDirectory . . . . .	734
27.275.2.5	FileIsSymlink . . . . .	734
27.275.2.6	FileSize . . . . .	734
27.275.2.7	FileTime . . . . .	734
27.275.2.8	FormatDateTime . . . . .	734
27.275.2.9	GetCurrentDateTime . . . . .	734
27.275.2.10	GetCurrentModuleFileName . . . . .	735
27.275.2.11	GetCurrentProcessFileName . . . . .	735
27.275.2.12	GetCurrentResourcesDirectory . . . . .	735
27.275.2.13	GetCurrentCWD . . . . .	735
27.275.2.14	GetHostName . . . . .	735
27.275.2.15	GetLastSystemError . . . . .	735
27.275.2.16	GetLocaleCharset . . . . .	735
27.275.2.17	GetPermissions . . . . .	735
27.275.2.18	GetTimezoneOffsetFromUTC . . . . .	735
27.275.2.19	MakeDirectory . . . . .	735
27.275.2.20	ParseDateTime . . . . .	736
27.275.2.21	ParseDateTime . . . . .	736
27.275.2.22	RemoveFile . . . . .	736
27.275.2.23	SetPermissions . . . . .	736
27.275.2.24	StrCaseCmp . . . . .	736
27.275.2.25	StrNCaseCmp . . . . .	736

27.275.2.28	StrSep	. 736
27.275.2.29	StrTokR	. 736
27.276	dcm::Table Class Reference	. 736
27.276.1	Detailed Description	. 737
27.276.2	Member Typedef Documentation	. 737
27.276.2.1	MapTableEntry	. 737
27.276.3	Constructor & Destructor Documentation	. 737
27.276.3.1	Table	. 737
27.276.3.2	~Table	. 737
27.276.4	Member Function Documentation	. 737
27.276.4.1	GetTableEntry	. 737
27.276.4.2	InsertEntry	. 737
27.276.5	Friends And Related Function Documentation	. 737
27.276.5.1	operator<<	. 737
27.277	dcm::TableEntry Class Reference	. 737
27.277.1	Detailed Description	. 738
27.277.2	Constructor & Destructor Documentation	. 738
27.277.2.1	TableEntry	. 738
27.277.2.2	~TableEntry	. 738
27.278	dcm::TableReader Class Reference	. 738
27.278.1	Detailed Description	. 739
27.278.2	Constructor & Destructor Documentation	. 739
27.278.2.1	TableReader	. 739
27.278.2.2	~TableReader	. 739
27.278.3	Member Function Documentation	. 739
27.278.3.1	CharacterDataHandler	. 739
27.278.3.2	EndElement	. 739
27.278.3.3	GetDefs	. 739
27.278.3.4	GetFilename	. 739
27.278.3.5	HandleIOD	. 739
27.278.3.6	HandleIODEntry	. 739
27.278.3.7	HandleMacro	. 739
27.278.3.8	HandleMacroEntry	. 739
27.278.3.9	HandleMacroEntryDescription	. 740
27.278.3.10	HandleModule	. 740
27.278.3.11	HandleModuleEntry	. 740
27.278.3.12	HandleModuleEntryDescription	. 740

27.278.3.1	<a href="#">HandleModuleInclude</a>	. 740
27.278.3.1	<a href="#">Read</a>	. 740
27.278.3.1	<a href="#">SetFilename</a>	. 740
27.278.3.1	<a href="#">StartElement</a>	. 740
27.279	<a href="#">dcm::network::TableRow Class Reference</a>	. 740
27.279	<a href="#">Constructor &amp; Destructor Documentation</a>	. 741
27.279.1	<a href="#">TableRow</a>	. 741
27.279.1.2	<a href="#">~TableRow</a>	. 741
27.279.2	<a href="#">Member Data Documentation</a>	. 741
27.279.2.1	<a href="#">transitions</a>	. 741
27.280	<a href="#">dcm::Tag Class Reference</a>	. 741
27.280	<a href="#">Detailed Description</a>	. 743
27.280	<a href="#">Constructor &amp; Destructor Documentation</a>	. 743
27.280.2	<a href="#">Tag</a>	. 743
27.280.2.2	<a href="#">Tag</a>	. 743
27.280.2.3	<a href="#">Tag</a>	. 744
27.280.3	<a href="#">Member Function Documentation</a>	. 744
27.280.3.1	<a href="#">GetElement</a>	. 744
27.280.3.2	<a href="#">GetElementTag</a>	. 744
27.280.3.3	<a href="#">GetGroup</a>	. 744
27.280.3.4	<a href="#">GetLength</a>	. 744
27.280.3.5	<a href="#">GetPrivateCreator</a>	. 744
27.280.3.6	<a href="#">IsGroupLength</a>	. 744
27.280.3.7	<a href="#">IsGroupXX</a>	. 745
27.280.3.8	<a href="#">IsIllegal</a>	. 745
27.280.3.9	<a href="#">IsPrivate</a>	. 745
27.280.3.10	<a href="#">IsPrivateCreator</a>	. 745
27.280.3.11	<a href="#">IsPublic</a>	. 745
27.280.3.12	<a href="#">operator!="</a>	. 745
27.280.3.13	<a href="#">operator&lt;</a>	. 745
27.280.3.14	<a href="#">operator&lt;=</a>	. 745
27.280.3.15	<a href="#">operator=</a>	. 745
27.280.3.16	<a href="#">operator==</a>	. 746
27.280.3.17	<a href="#">operator[]</a>	. 746
27.280.3.18	<a href="#">operator[]</a>	. 746
27.280.3.19	<a href="#">PrintAsContinuousString</a>	. 746
27.280.3.20	<a href="#">PrintAsContinuousUpperCaseString</a>	. 746

27.280.3.21	PrintAsPipeSeparatedString	746
27.280.3.22	Read	746
27.280.3.23	ReadFromCommaSeparatedString	746
27.280.3.24	ReadFromContinuousString	746
27.280.3.25	ReadFromPipeSeparatedString	747
27.280.3.26	SetElement	747
27.280.3.27	SetElementTag	747
27.280.3.28	SetElementTag	747
27.280.3.29	SetGroup	747
27.280.3.30	SetPrivateCreator	747
27.280.3.31	Write	747
27.280.4	Friends And Related Function Documentation	748
27.280.4.1	operator<<	748
27.280.4.2	operator>>	748
27.280.5	Member Data Documentation	748
27.280.5.1	bytes	748
27.280.5.2	tag	748
27.280.5.3	tags	748
27.281	gdcmm::TagPath Class Reference	748
27.281.1	Detailed Description	748
27.281.2	Constructor & Destructor Documentation	749
27.281.2.1	TagPath	749
27.281.2.2	~TagPath	749
27.281.3	Member Function Documentation	749
27.281.3.1	ConstructFromString	749
27.281.3.2	ConstructFromTagList	749
27.281.3.3	IsValid	749
27.281.3.4	Print	749
27.281.3.5	Push	749
27.281.3.6	Push	749
27.282	gdcmm::Testing Class Reference	749
27.282.1	Detailed Description	750
27.282.2	Member Typedef Documentation	750
27.282.2.1	IMD5DataImagesType	750
27.282.2.2	MediaStorageDataFilesType	751
27.282.3	Constructor & Destructor Documentation	751
27.282.3.1	Testing	751

27.282.3.2~Testing	. 751
27.282.4Member Function Documentation	. 751
27.282.4.1ComputeFileMD5	. 751
27.282.4.2ComputeMD5	. 751
27.282.4.3GetDataExtraRoot	. 751
27.282.4.4GetDataRoot	. 751
27.282.4.5GetFileName	. 751
27.282.4.6GetFileNames	. 752
27.282.4.7GetLossyFlagFromFile	. 752
27.282.4.8GetMD5DataImage	. 752
27.282.4.9GetMD5DataImages	. 752
27.282.4.10GetMD5FromBrokenFile	. 752
27.282.4.11GetMD5FromFile	. 752
27.282.4.12GetMediaStorageDataFile	. 752
27.282.4.13GetMediaStorageDataFiles	. 752
27.282.4.14GetMediaStorageFromFile	. 752
27.282.4.15GetNumberOfFileNames	. 752
27.282.4.16GetNumberOfMD5DataImages	. 752
27.282.4.17GetNumberOfMediaStorageDataFiles	. 752
27.282.4.18GetPixelSpacingDataRoot	. 752
27.282.4.19GetSelectedPrivateGroupOffsetFromFile	. 753
27.282.4.20GetSelectedTagsOffsetFromFile	. 753
27.282.4.21GetSourceDirectory	. 753
27.282.4.22GetStreamOffsetFromFile	. 753
27.282.4.23GetTempDirectory	. 753
27.282.4.24GetTempDirectoryW	. 753
27.282.4.25GetTempFilename	. 753
27.282.4.26GetTempFilenameW	. 753
27.282.4.27Print	. 754
27.283dcm::Trace Class Reference	. 754
27.283.1Detailed Description	. 755
27.283.2Constructor & Destructor Documentation	. 755
27.283.2.1Trace	. 755
27.283.2.2~Trace	. 755
27.283.3Member Function Documentation	. 755
27.283.3.1DebugOff	. 755
27.283.3.2DebugOn	. 755

27.283.3.3ErrorOff . . . . .	755
27.283.3.4ErrorOn . . . . .	755
27.283.3.5GetDebugFlag . . . . .	755
27.283.3.6GetDebugStream . . . . .	755
27.283.3.7GetErrorFlag . . . . .	755
27.283.3.8GetErrorStream . . . . .	756
27.283.3.9GetStream . . . . .	756
27.283.3.10GetWarningFlag . . . . .	756
27.283.3.10GetWarningStream . . . . .	756
27.283.3.12SetDebug . . . . .	756
27.283.3.13SetDebugStream . . . . .	756
27.283.3.13SetError . . . . .	756
27.283.3.13SetErrorStream . . . . .	756
27.283.3.13SetStream . . . . .	756
27.283.3.13SetStreamToFile . . . . .	756
27.283.3.13SetWarning . . . . .	756
27.283.3.13SetWarningStream . . . . .	757
27.283.3.20WarningOff . . . . .	757
27.283.3.20WarningOn . . . . .	757
27.284dcm::TransferSyntax Class Reference . . . . .	757
27.284.1Detailed Description . . . . .	759
27.284.2Member Enumeration Documentation . . . . .	759
27.284.2.1NegociatedType . . . . .	759
27.284.2.2TSType . . . . .	759
27.284.3Constructor & Destructor Documentation . . . . .	760
27.284.3.1TransferSyntax . . . . .	760
27.284.4Member Function Documentation . . . . .	760
27.284.4.1CanStoreLossy . . . . .	760
27.284.4.2GetNegociatedType . . . . .	760
27.284.4.3GetString . . . . .	760
27.284.4.4GetSwapCode . . . . .	760
27.284.4.5GetTSString . . . . .	760
27.284.4.6GetTSType . . . . .	760
27.284.4.7IsEncapsulated . . . . .	760
27.284.4.8IsEncoded . . . . .	761
27.284.4.9IsExplicit . . . . .	761
27.284.4.10IsImplicit . . . . .	761

27.284.4.11\$Lossless	. 761
27.284.4.11\$Lossy	. 761
27.284.4.11\$Valid	. 761
27.284.4.1operator TSType	. 761
27.284.5Friends And Related Function Documentation	. 761
27.284.5.1operator<<	. 761
27.285gdcmm::network::TransferSyntaxSub Class Reference	. 761
27.285.1Detailed Description	. 761
27.285.2Constructor & Destructor Documentation	. 762
27.285.2.1TransferSyntaxSub	. 762
27.285.3Member Function Documentation	. 762
27.285.3.1GetName	. 762
27.285.3.2operator==	. 762
27.285.3.3Print	. 762
27.285.3.4Read	. 762
27.285.3.5SetName	. 762
27.285.3.6SetNameFromUID	. 762
27.285.3.7Size	. 762
27.285.3.8Write	. 762
27.286gdcmm::network::Transition Struct Reference	. 762
27.286.1Constructor & Destructor Documentation	. 763
27.286.1.1Transition	. 763
27.286.1.2~Transition	. 763
27.286.1.3Transition	. 763
27.286.2Member Function Documentation	. 763
27.286.2.1MakeNew	. 763
27.286.3Member Data Documentation	. 763
27.286.3.1mAction	. 763
27.286.3.2mEnd	. 763
27.287gdcmm::Type Class Reference	. 764
27.287.1Detailed Description	. 764
27.287.2Member Enumeration Documentation	. 765
27.287.2.1TypeType	. 765
27.287.3Constructor & Destructor Documentation	. 765
27.287.3.1Type	. 765
27.287.4Member Function Documentation	. 765
27.287.4.1GetTypeString	. 765



27.287.4.2GetTypeType . . . . .	765
27.287.4.3operator TypeType . . . . .	765
27.287.5Friends And Related Function Documentation . . . . .	765
27.287.5.1operator<< . . . . .	765
27.288dcm::UI Struct Reference . . . . .	765
27.288.1Friends And Related Function Documentation . . . . .	766
27.288.1.1operator<< . . . . .	766
27.288.2Member Data Documentation . . . . .	766
27.288.2.1Internal . . . . .	766
27.289dcm::UIDGenerator Class Reference . . . . .	766
27.289.1Detailed Description . . . . .	766
27.289.2Constructor & Destructor Documentation . . . . .	767
27.289.2.1UIDGenerator . . . . .	767
27.289.3Member Function Documentation . . . . .	767
27.289.3.1Generate . . . . .	767
27.289.3.2GenerateUUID . . . . .	767
27.289.3.3GetGDCMUID . . . . .	767
27.289.3.4GetRoot . . . . .	767
27.289.3.5IsValid . . . . .	767
27.289.3.6SetRoot . . . . .	768
27.290dcm::UIDs Class Reference . . . . .	768
27.290.1Detailed Description . . . . .	772
27.290.2Member Typedef Documentation . . . . .	772
27.290.2.1TransferSyntaxStringsType . . . . .	772
27.290.3Member Enumeration Documentation . . . . .	772
27.290.3.1TSName . . . . .	772
27.290.3.2TSType . . . . .	779
27.290.4Member Function Documentation . . . . .	785
27.290.4.1GetName . . . . .	785
27.290.4.2GetNumberOfTransferSyntaxStrings . . . . .	786
27.290.4.3GetString . . . . .	786
27.290.4.4GetTransferSyntaxString . . . . .	786
27.290.4.5GetTransferSyntaxStrings . . . . .	786
27.290.4.6GetUIDName . . . . .	786
27.290.4.7GetUIDString . . . . .	786
27.290.4.8operator TSType . . . . .	786
27.290.4.9SetFromUID . . . . .	786

27.294	dcm::network::ULAction Class Reference	. 786
27.291.1	Detailed Description	. 788
27.291.2	Constructor & Destructor Documentation	. 788
27.291.2.1	ULAction	. 788
27.291.2.2	~ULAction	. 788
27.291.3	Member Function Documentation	. 788
27.291.3.1	PerformAction	. 788
27.294	dcm::network::ULActionAA1 Class Reference	. 789
27.292.1	Member Function Documentation	. 789
27.292.1.1	PerformAction	. 789
27.294	dcm::network::ULActionAA2 Class Reference	. 790
27.293.1	Member Function Documentation	. 790
27.293.1.1	PerformAction	. 790
27.294	dcm::network::ULActionAA3 Class Reference	. 791
27.294.1	Member Function Documentation	. 791
27.294.1.1	PerformAction	. 791
27.294	dcm::network::ULActionAA4 Class Reference	. 792
27.295.1	Member Function Documentation	. 792
27.295.1.1	PerformAction	. 792
27.294	dcm::network::ULActionAA5 Class Reference	. 793
27.296.1	Member Function Documentation	. 793
27.296.1.1	PerformAction	. 793
27.294	dcm::network::ULActionAA6 Class Reference	. 794
27.297.1	Member Function Documentation	. 794
27.297.1.1	PerformAction	. 794
27.294	dcm::network::ULActionAA7 Class Reference	. 795
27.298.1	Member Function Documentation	. 795
27.298.1.1	PerformAction	. 795
27.294	dcm::network::ULActionAA8 Class Reference	. 796
27.299.1	Member Function Documentation	. 796
27.299.1.1	PerformAction	. 796
27.300	dcm::network::ULActionAE1 Class Reference	. 797
27.300.1	Member Function Documentation	. 797
27.300.1.1	PerformAction	. 797
27.304	dcm::network::ULActionAE2 Class Reference	. 798
27.301.1	Member Function Documentation	. 798
27.301.1.1	PerformAction	. 798

27.302dcm::network::ULActionAE3 Class Reference . . . . .	799
27.302.1Member Function Documentation . . . . .	799
27.302.1.1PerformAction . . . . .	799
27.303dcm::network::ULActionAE4 Class Reference . . . . .	800
27.303.1Member Function Documentation . . . . .	800
27.303.1.1PerformAction . . . . .	800
27.304dcm::network::ULActionAE5 Class Reference . . . . .	801
27.304.1Member Function Documentation . . . . .	801
27.304.1.1PerformAction . . . . .	801
27.305dcm::network::ULActionAE6 Class Reference . . . . .	802
27.305.1Member Function Documentation . . . . .	802
27.305.1.1PerformAction . . . . .	802
27.306dcm::network::ULActionAE7 Class Reference . . . . .	803
27.306.1Member Function Documentation . . . . .	803
27.306.1.1PerformAction . . . . .	803
27.307dcm::network::ULActionAE8 Class Reference . . . . .	804
27.307.1Member Function Documentation . . . . .	804
27.307.1.1PerformAction . . . . .	804
27.308dcm::network::ULActionAR1 Class Reference . . . . .	805
27.308.1Member Function Documentation . . . . .	805
27.308.1.1PerformAction . . . . .	805
27.309dcm::network::ULActionAR10 Class Reference . . . . .	806
27.309.1Member Function Documentation . . . . .	806
27.309.1.1PerformAction . . . . .	806
27.310dcm::network::ULActionAR2 Class Reference . . . . .	807
27.310.1Member Function Documentation . . . . .	807
27.310.1.1PerformAction . . . . .	807
27.311dcm::network::ULActionAR3 Class Reference . . . . .	808
27.311.1Member Function Documentation . . . . .	808
27.311.1.1PerformAction . . . . .	808
27.312dcm::network::ULActionAR4 Class Reference . . . . .	809
27.312.1Member Function Documentation . . . . .	809
27.312.1.1PerformAction . . . . .	809
27.313dcm::network::ULActionAR5 Class Reference . . . . .	810
27.313.1Member Function Documentation . . . . .	810
27.313.1.1PerformAction . . . . .	810
27.314dcm::network::ULActionAR6 Class Reference . . . . .	811

27.314.1Member Function Documentation . . . . .	811
27.314.1.1PerformAction . . . . .	811
27.315dcm::network::ULActionAR7 Class Reference . . . . .	812
27.315.1Member Function Documentation . . . . .	812
27.315.1.1PerformAction . . . . .	812
27.316dcm::network::ULActionAR8 Class Reference . . . . .	813
27.316.1Member Function Documentation . . . . .	813
27.316.1.1PerformAction . . . . .	813
27.317dcm::network::ULActionAR9 Class Reference . . . . .	814
27.317.1Member Function Documentation . . . . .	814
27.317.1.1PerformAction . . . . .	814
27.318dcm::network::ULActionDT1 Class Reference . . . . .	815
27.318.1Member Function Documentation . . . . .	815
27.318.1.1PerformAction . . . . .	815
27.319dcm::network::ULActionDT2 Class Reference . . . . .	816
27.319.1Member Function Documentation . . . . .	816
27.319.1.1PerformAction . . . . .	816
27.320dcm::network::ULBasicCallback Class Reference . . . . .	817
27.320.1Detailed Description . . . . .	818
27.320.2Constructor & Destructor Documentation . . . . .	818
27.320.2.1ULBasicCallback . . . . .	818
27.320.2.2~ULBasicCallback . . . . .	818
27.320.3Member Function Documentation . . . . .	818
27.320.3.1GetDataSets . . . . .	818
27.320.3.2GetResponses . . . . .	818
27.320.3.3HandleDataSet . . . . .	818
27.320.3.4HandleResponse . . . . .	818
27.321dcm::network::ULConnection Class Reference . . . . .	818
27.321.1Detailed Description . . . . .	819
27.321.2Constructor & Destructor Documentation . . . . .	819
27.321.2.1ULConnection . . . . .	819
27.321.2.2~ULConnection . . . . .	819
27.321.3Member Function Documentation . . . . .	819
27.321.3.1AddAcceptedPresentationContext . . . . .	819
27.321.3.2FindContext . . . . .	820
27.321.3.3GetAcceptedPresentationContexts . . . . .	820
27.321.3.4GetAcceptedPresentationContexts . . . . .	820

27.321.3.5	GetConnectionInfo	820
27.321.3.6	GetMaxPDUSize	820
27.321.3.7	GetPresentationContextACByID	820
27.321.3.8	GetPresentationContextIDFromPresentationContext	820
27.321.3.9	GetPresentationContextRQByID	820
27.321.3.10	GetPresentationContexts	820
27.321.3.10	GetProtocol	820
27.321.3.10	GetState	820
27.321.3.10	GetTimer	820
27.321.3.11	InitializeConnection	820
27.321.3.11	InitializeIncomingConnection	820
27.321.3.11	SetMaxPDUSize	820
27.321.3.11	SetPresentationContexts	820
27.321.3.11	SetPresentationContexts	820
27.321.3.11	SetState	820
27.321.3.20	StopProtocol	820
27.321.4	Friends And Related Function Documentation	820
27.321.4.1	ULActionAE6	821
27.321.4.2	ULConnectionManager	821
27.322	dcm::network::ULConnectionCallback Class Reference	821
27.322.1	Detailed Description	822
27.322.2	Constructor & Destructor Documentation	822
27.322.2.1	ULConnectionCallback	822
27.322.2.2	~ULConnectionCallback	822
27.322.3	Member Function Documentation	822
27.322.3.1	DataSetHandled	822
27.322.3.2	DataSetHandles	822
27.322.3.3	HandleDataSet	822
27.322.3.4	HandleResponse	822
27.322.3.5	ResetHandledDataSet	822
27.322.3.6	SetImplicitFlag	822
27.322.4	Member Data Documentation	822
27.322.4.1	Implicit	822
27.323	dcm::network::ULConnectionInfo Class Reference	822
27.323.1	Detailed Description	823
27.323.2	Constructor & Destructor Documentation	823
27.323.2.1	ULConnectionInfo	823

27.323.3	Member Function Documentation	823
27.323.3.1	GetCalledAETitle	823
27.323.3.2	GetCalledComputerName	823
27.323.3.3	GetCalledIPAddress	823
27.323.3.4	GetCalledIPPort	823
27.323.3.5	GetCallingAETitle	823
27.323.3.6	GetMaxPDULength	823
27.323.3.7	Initialize	823
27.323.3.8	SetMaxPDULength	823
27.324	dcm::network::ULConnectionManager Class Reference	824
27.324.1	Detailed Description	825
27.324.2	Constructor & Destructor Documentation	825
27.324.2.1	ULConnectionManager	825
27.324.2.2	~ULConnectionManager	825
27.324.3	Member Function Documentation	825
27.324.3.1	BreakConnection	825
27.324.3.2	BreakConnectionNow	825
27.324.3.3	EstablishConnection	825
27.324.3.4	EstablishConnectionMove	826
27.324.3.5	SendEcho	826
27.324.3.6	SendFind	826
27.324.3.7	SendFind	826
27.324.3.8	SendMove	826
27.324.3.9	SendMove	826
27.324.3.10	SendStore	826
27.324.3.11	SendStore	826
27.325	dcm::network::ULEvent Class Reference	826
27.325.1	Detailed Description	826
27.325.2	Constructor & Destructor Documentation	827
27.325.2.1	ULEvent	827
27.325.2.2	ULEvent	827
27.325.2.3	~ULEvent	827
27.325.3	Member Function Documentation	827
27.325.3.1	GetEvent	827
27.325.3.2	GetPDUs	827
27.325.3.3	SetEvent	827
27.325.3.4	SetPDU	827

27.326	dcm::network::ULTransitionTable Class Reference	827
27.326.1	Detailed Description	827
27.326.2	Constructor & Destructor Documentation	828
27.326.2.1	ULTransitionTable	828
27.326.3	Member Function Documentation	828
27.326.3.1	HandleEvent	828
27.326.3.2	PrintTable	828
27.327	dcm::network::ULWritingCallback Class Reference	828
27.327.1	Constructor & Destructor Documentation	829
27.327.1.1	ULWritingCallback	829
27.327.1.2	~ULWritingCallback	829
27.327.2	Member Function Documentation	829
27.327.2.1	HandleDataSet	829
27.327.2.2	HandleResponse	829
27.327.2.3	SetDirectory	829
27.328	dcm::UNExplicitDataElement Class Reference	830
27.328.1	Detailed Description	831
27.328.2	Member Function Documentation	831
27.328.2.1	GetLength	831
27.328.2.2	Read	831
27.328.2.3	ReadPreValue	831
27.328.2.4	ReadValue	831
27.328.2.5	ReadWithLength	831
27.329	dcm::UNExplicitImplicitDataElement Class Reference	831
27.329.1	Detailed Description	833
27.329.2	Member Function Documentation	833
27.329.2.1	GetLength	833
27.329.2.2	Read	833
27.329.2.3	ReadPreValue	833
27.329.2.4	ReadValue	833
27.330	dcm::Unpacker12Bits Class Reference	833
27.330.1	Detailed Description	833
27.330.2	Member Function Documentation	834
27.330.2.1	Pack	834
27.330.2.2	Unpack	834
27.331	dcm::Usage Class Reference	834
27.331.1	Detailed Description	835

27.331.2	Member Enumeration Documentation	835
27.331.2.1	UsageType	835
27.331.3	Constructor & Destructor Documentation	835
27.331.3.1	Usage	835
27.331.4	Member Function Documentation	835
27.331.4.1	GetUsageString	835
27.331.4.2	GetUsageType	835
27.331.4.3	operator UsageType	835
27.331.5	Friends And Related Function Documentation	835
27.331.5.1	operator<<	836
27.332	dcm::UserEvent Class Reference	836
27.333	dcm::network::UserInformation Class Reference	837
27.333.1	Detailed Description	837
27.333.2	Constructor & Destructor Documentation	837
27.333.2.1	UserInformation	837
27.333.2.2	~UserInformation	837
27.333.3	Member Function Documentation	837
27.333.3.1	AddRoleSelectionSub	837
27.333.3.2	AddSOPClassExtendedNegociationSub	837
27.333.3.3	GetMaximumLengthSub	837
27.333.3.4	GetMaximumLengthSub	838
27.333.3.5	operator=	838
27.333.3.6	Print	838
27.333.3.7	Read	838
27.333.3.8	Size	838
27.333.3.9	Write	838
27.334	dcm::UUIDGenerator Class Reference	838
27.334.1	Detailed Description	838
27.334.2	Member Function Documentation	838
27.334.2.1	Generate	838
27.334.2.2	IsValid	838
27.335	dcm::Validate Class Reference	839
27.335.1	Detailed Description	839
27.335.2	Constructor & Destructor Documentation	839
27.335.2.1	Validate	839
27.335.2.2	~Validate	840
27.335.3	Member Function Documentation	840



27.335.3.1GetValidatedFile . . . . .	840
27.335.3.2SetFile . . . . .	840
27.335.3.3Validation . . . . .	840
27.335.4Member Data Documentation . . . . .	840
27.335.4.1F . . . . .	840
27.335.4.2V . . . . .	840
27.336dcm::Value Class Reference . . . . .	840
27.336.1Detailed Description . . . . .	841
27.336.2Constructor & Destructor Documentation . . . . .	841
27.336.2.1Value . . . . .	841
27.336.2.2~Value . . . . .	841
27.336.3Member Function Documentation . . . . .	841
27.336.3.1Clear . . . . .	842
27.336.3.2GetLength . . . . .	842
27.336.3.3operator== . . . . .	842
27.336.3.4SetLength . . . . .	842
27.336.3.5SetLengthOnly . . . . .	842
27.336.4Friends And Related Function Documentation . . . . .	842
27.336.4.1DataElement . . . . .	842
27.337dcm::ValueIO< TDE, TSwap, TType > Class Template Reference . . . . .	842
27.337.1Detailed Description . . . . .	842
27.337.2Member Function Documentation . . . . .	843
27.337.2.1Read . . . . .	843
27.337.2.2Write . . . . .	843
27.338dcm::Version Class Reference . . . . .	843
27.338.1Detailed Description . . . . .	843
27.338.2Constructor & Destructor Documentation . . . . .	843
27.338.2.1Version . . . . .	843
27.338.2.2~Version . . . . .	843
27.338.3Member Function Documentation . . . . .	843
27.338.3.1GetBuildVersion . . . . .	844
27.338.3.2GetMajorVersion . . . . .	844
27.338.3.3GetMinorVersion . . . . .	844
27.338.3.4GetVersion . . . . .	844
27.338.3.5Print . . . . .	844
27.338.4Friends And Related Function Documentation . . . . .	844
27.338.4.1operator<< . . . . .	844

27.339	dcm::VL Class Reference . . . . .	844
27.339.1	Detailed Description . . . . .	845
27.339.2	Member Typedef Documentation . . . . .	845
27.339.2.1	Type . . . . .	845
27.339.3	Constructor & Destructor Documentation . . . . .	845
27.339.3.1	VL . . . . .	845
27.339.4	Member Function Documentation . . . . .	845
27.339.4.1	GetLength . . . . .	845
27.339.4.2	GetVL16Max . . . . .	845
27.339.4.3	GetVL32Max . . . . .	845
27.339.4.4	IsOdd . . . . .	845
27.339.4.5	IsUndefined . . . . .	846
27.339.4.6	operator uint32_t . . . . .	846
27.339.4.7	operator++ . . . . .	846
27.339.4.8	operator++ . . . . .	846
27.339.4.9	operator+= . . . . .	846
27.339.4.10	Read . . . . .	846
27.339.4.11	Read16 . . . . .	846
27.339.4.12	SetToUndefined . . . . .	846
27.339.4.13	Write . . . . .	846
27.339.4.14	Write16 . . . . .	846
27.339.5	Friends And Related Function Documentation . . . . .	846
27.339.5.1	operator<< . . . . .	846
27.340	dcm::VM Class Reference . . . . .	846
27.340.1	Detailed Description . . . . .	848
27.340.2	Member Enumeration Documentation . . . . .	848
27.340.2.1	VMType . . . . .	848
27.340.3	Constructor & Destructor Documentation . . . . .	849
27.340.3.1	VM . . . . .	849
27.340.4	Member Function Documentation . . . . .	849
27.340.4.1	Compatible . . . . .	849
27.340.4.2	GetIndex . . . . .	849
27.340.4.3	GetLength . . . . .	849
27.340.4.4	GetNumberOfElementsFromArray . . . . .	849
27.340.4.5	GetVMString . . . . .	849
27.340.4.6	GetVMType . . . . .	850
27.340.4.7	GetVMTypeFromLength . . . . .	850

27.340.4.8IsValid	850
27.340.4.9operator VMType	850
27.340.5Friends And Related Function Documentation	850
27.340.5.1operator<<	850
27.340dcm::VMToLength< T > Struct Template Reference	850
27.340dcm::VR Class Reference	850
27.342.1Detailed Description	852
27.342.2Member Enumeration Documentation	852
27.342.2.1VRType	852
27.342.3Constructor & Destructor Documentation	854
27.342.3.1VR	854
27.342.4Member Function Documentation	854
27.342.4.1CanDisplay	854
27.342.4.2Compatible	854
27.342.4.3GetLength	854
27.342.4.4GetLength	854
27.342.4.5GetSize	854
27.342.4.6GetSizeof	854
27.342.4.7GetVRString	854
27.342.4.8GetVRStringFromFile	854
27.342.4.9GetVRType	854
27.342.4.10GetVRTypeFromFile	854
27.342.4.11ASCII	854
27.342.4.12ASCII2	854
27.342.4.13Binary	854
27.342.4.14Binary2	854
27.342.4.15Dual	854
27.342.4.16Swap	854
27.342.4.17Valid	854
27.342.4.18Valid	854
27.342.4.19VRFile	854
27.342.4.20operator VRType	855
27.342.4.21Read	855
27.342.4.22Write	855
27.342.5Friends And Related Function Documentation	855
27.342.5.1operator<<	855
27.340dcm::VR16ExplicitDataElement Class Reference	855

27.343.1Detailed Description . . . . .	856
27.343.2Member Function Documentation . . . . .	856
27.343.2.1GetLength . . . . .	856
27.343.2.2Read . . . . .	856
27.343.2.3ReadPreValue . . . . .	857
27.343.2.4ReadValue . . . . .	857
27.343.2.5ReadWithLength . . . . .	857
27.344dcm::VRToEncoding< T > Struct Template Reference . . . . .	857
27.345dcm::VRToType< T > Struct Template Reference . . . . .	857
27.345.1Detailed Description . . . . .	857
27.346dcm::VRVLSize< T > Class Template Reference . . . . .	858
27.347dcm::VRVLSize< 0 > Class Template Reference . . . . .	858
27.347.1Member Function Documentation . . . . .	858
27.347.1.1Read . . . . .	858
27.347.1.2Write . . . . .	858
27.348dcm::VRVLSize< 1 > Class Template Reference . . . . .	858
27.348.1Member Function Documentation . . . . .	858
27.348.1.1Read . . . . .	858
27.348.1.2Write . . . . .	858
27.349tkGDCMImageReader Class Reference . . . . .	859
27.349.1Detailed Description . . . . .	861
27.349.2Constructor & Destructor Documentation . . . . .	861
27.349.2.1tkGDCMImageReader . . . . .	861
27.349.2.2~tkGDCMImageReader . . . . .	862
27.349.3Member Function Documentation . . . . .	862
27.349.3.1CanReadFile . . . . .	862
27.349.3.2ExecuteData . . . . .	862
27.349.3.3ExecuteInformation . . . . .	862
27.349.3.4FillMedicalImageInformation . . . . .	862
27.349.3.5GetDescriptiveName . . . . .	862
27.349.3.6GetFileExtensions . . . . .	862
27.349.3.7GetIconImage . . . . .	862
27.349.3.8GetOverlay . . . . .	862
27.349.3.9LoadSingleFile . . . . .	862
27.349.3.10New . . . . .	862
27.349.3.11PrintSelf . . . . .	862
27.349.3.12requestDataCompat . . . . .	862

27.349.3.18	RequestInformationCompat	862
27.349.3.19	SetCurve	862
27.349.3.19	SetFileNames	863
27.349.3.19	SetFilePattern	863
27.349.3.19	SetFilePrefix	863
27.349.3.19	SetMedicalImageProperties	863
27.349.3.19	tkBooleanMacro	863
27.349.3.20	tkBooleanMacro	863
27.349.3.21	tkBooleanMacro	863
27.349.3.22	tkBooleanMacro	863
27.349.3.23	tkBooleanMacro	863
27.349.3.24	tkGetMacro	863
27.349.3.25	tkGetMacro	863
27.349.3.26	tkGetMacro	863
27.349.3.27	tkGetMacro	863
27.349.3.28	tkGetMacro	863
27.349.3.29	tkGetMacro	863
27.349.3.30	tkGetMacro	863
27.349.3.31	tkGetMacro	863
27.349.3.32	tkGetMacro	863
27.349.3.33	tkGetMacro	863
27.349.3.34	tkGetMacro	863
27.349.3.35	tkGetObjectMacro	863
27.349.3.36	tkGetObjectMacro	863
27.349.3.37	tkGetObjectMacro	863
27.349.3.38	tkGetObjectMacro	864
27.349.3.39	tkGetStringMacro	864
27.349.3.40	tkGetStringMacro	864
27.349.3.41	tkGetVector3Macro	864
27.349.3.42	tkGetVector6Macro	864
27.349.3.43	tkSetMacro	864
27.349.3.44	tkSetMacro	864
27.349.3.45	tkSetMacro	864
27.349.3.46	tkSetMacro	864
27.349.3.47	tkSetVector6Macro	864
27.349.3.48	tkTypeRevisionMacro	864
27.349.4	Member Data Documentation	864

27.349.4.1	ApplyInverseVideo	864
27.349.4.2	ApplyLookupTable	864
27.349.4.3	ApplyPlanarConfiguration	864
27.349.4.4	ApplyShiftScale	864
27.349.4.5	ApplyYBRToRGB	864
27.349.4.6	Curve	864
27.349.4.7	DirectionCosines	864
27.349.4.8	FileNames	864
27.349.4.9	ForceRescale	864
27.349.4.10	IconDataScalarType	864
27.349.4.11	IconImageDataExtent	864
27.349.4.12	IconNumberOfScalarComponents	864
27.349.4.13	IconImageFormat	864
27.349.4.14	IconImageOrientationPatient	864
27.349.4.15	IconImagePositionPatient	865
27.349.4.16	LoadIconImage	865
27.349.4.17	LoadOverlays	865
27.349.4.18	LossyFlag	865
27.349.4.19	MedicalImageProperties	865
27.349.4.20	NumberOfIconImages	865
27.349.4.21	NumberOfOverlays	865
27.349.4.22	PlanarConfiguration	865
27.349.4.23	Scale	865
27.349.4.24	Shift	865
27.350	tkGDCMImageReader2 Class Reference	865
27.350.1	Detailed Description	868
27.350.2	Constructor & Destructor Documentation	868
27.350.2.1	tkGDCMImageReader2	868
27.350.2.2	~tkGDCMImageReader2	868
27.350.3	Member Function Documentation	868
27.350.3.1	CanReadFile	868
27.350.3.2	FillMedicalImageInformation	868
27.350.3.3	GetDescriptiveName	868
27.350.3.4	GetFileExtensions	868
27.350.3.5	GetIconImage	868
27.350.3.6	GetIconImagePort	868
27.350.3.7	GetOverlay	868

27.350.3.8GetOverlayPort . . . . .	868
27.350.3.9LoadSingleFile . . . . .	868
27.350.3.10New . . . . .	868
27.350.3.11PrintSelf . . . . .	868
27.350.3.12ProcessRequest . . . . .	868
27.350.3.13RequestData . . . . .	868
27.350.3.14RequestDataCompat . . . . .	868
27.350.3.15RequestInformation . . . . .	868
27.350.3.16RequestInformationCompat . . . . .	869
27.350.3.17SetCurve . . . . .	869
27.350.3.18SetFilePattern . . . . .	869
27.350.3.19SetFilePrefix . . . . .	869
27.350.3.20SetMedicalImageProperties . . . . .	869
27.350.3.21tkBooleanMacro . . . . .	869
27.350.3.22tkBooleanMacro . . . . .	869
27.350.3.23tkBooleanMacro . . . . .	869
27.350.3.24tkBooleanMacro . . . . .	869
27.350.3.25tkBooleanMacro . . . . .	869
27.350.3.26tkGetMacro . . . . .	869
27.350.3.27tkGetMacro . . . . .	869
27.350.3.28tkGetMacro . . . . .	869
27.350.3.29tkGetMacro . . . . .	869
27.350.3.30tkGetMacro . . . . .	869
27.350.3.31tkGetMacro . . . . .	869
27.350.3.32tkGetMacro . . . . .	869
27.350.3.33tkGetMacro . . . . .	869
27.350.3.34tkGetMacro . . . . .	869
27.350.3.35tkGetMacro . . . . .	869
27.350.3.36tkGetMacro . . . . .	869
27.350.3.37tkGetObjectMacro . . . . .	869
27.350.3.38tkGetObjectMacro . . . . .	869
27.350.3.39tkGetStringMacro . . . . .	869
27.350.3.40tkGetStringMacro . . . . .	869
27.350.3.41tkGetVector3Macro . . . . .	869
27.350.3.42tkGetVector6Macro . . . . .	870
27.350.3.43tkSetMacro . . . . .	870
27.350.3.44tkSetMacro . . . . .	870

27.350.3.46kSetMacro . . . . .	870
27.350.3.46kSetMacro . . . . .	870
27.350.3.47kSetVector6Macro . . . . .	870
27.350.3.48kTypeRevisionMacro . . . . .	870
27.350.4Member Data Documentation . . . . .	870
27.350.4.1ApplyInverseVideo . . . . .	870
27.350.4.2ApplyLookupTable . . . . .	870
27.350.4.3ApplyPlanarConfiguration . . . . .	870
27.350.4.4ApplyShiftScale . . . . .	870
27.350.4.5ApplyYBRToRGB . . . . .	870
27.350.4.6Curve . . . . .	870
27.350.4.7DirectionCosines . . . . .	870
27.350.4.8ForceRescale . . . . .	870
27.350.4.9IconDataScalarType . . . . .	870
27.350.4.10IconImageDataExtent . . . . .	870
27.350.4.11IconNumberOfScalarComponents . . . . .	870
27.350.4.12ImageFormat . . . . .	870
27.350.4.13ImageOrientationPatient . . . . .	870
27.350.4.14ImagePositionPatient . . . . .	870
27.350.4.15LoadIconImage . . . . .	870
27.350.4.16LoadOverlays . . . . .	870
27.350.4.17LossyFlag . . . . .	870
27.350.4.18NumberOfIconImages . . . . .	870
27.350.4.19NumberOfOverlays . . . . .	871
27.350.4.20PlanarConfiguration . . . . .	871
27.350.4.21Scale . . . . .	871
27.350.4.22Shift . . . . .	871
27.351.1vtkGDCMImageWriter Class Reference . . . . .	871
27.351.1Detailed Description . . . . .	873
27.351.2Member Enumeration Documentation . . . . .	873
27.351.2.1CompressionTypes . . . . .	873
27.351.3Constructor & Destructor Documentation . . . . .	873
27.351.3.1vtkGDCMImageWriter . . . . .	873
27.351.3.2~vtkGDCMImageWriter . . . . .	873
27.351.4Member Function Documentation . . . . .	873
27.351.4.1GetDescriptiveName . . . . .	873
27.351.4.2GetFileExtensions . . . . .	873



27.351.4.3GetFileName	873
27.351.4.4New	873
27.351.4.5PrintSelf	874
27.351.4.6SetDirectionCosines	874
27.351.4.7SetDirectionCosinesFromImageOrientationPatient	874
27.351.4.8SetFileNames	874
27.351.4.9SetMedicalImageProperties	874
27.351.4.10tkBooleanMacro	874
27.351.4.11tkBooleanMacro	874
27.351.4.12tkGetMacro	874
27.351.4.13tkGetMacro	874
27.351.4.14tkGetMacro	874
27.351.4.15tkGetMacro	874
27.351.4.16tkGetMacro	874
27.351.4.17tkGetMacro	874
27.351.4.18tkGetMacro	874
27.351.4.19tkGetObjectMacro	874
27.351.4.20tkGetObjectMacro	874
27.351.4.21tkGetObjectMacro	874
27.351.4.22tkGetStringMacro	874
27.351.4.23tkGetStringMacro	875
27.351.4.24tkSetMacro	875
27.351.4.25tkSetMacro	875
27.351.4.26tkSetMacro	875
27.351.4.27tkSetMacro	875
27.351.4.28tkSetMacro	875
27.351.4.29tkSetMacro	875
27.351.4.30tkSetMacro	875
27.351.4.31tkSetStringMacro	875
27.351.4.32tkSetStringMacro	875
27.351.4.33tkTypeRevisionMacro	875
27.351.4.34Write	875
27.351.4.35WriteGDCMData	875
27.351.4.36WriteSlice	875
27.352tkGDCMMedicalImageProperties Class Reference	875
27.352.1Constructor & Destructor Documentation	877
27.352.1.1tkGDCMMedicalImageProperties	877

27.352.1.2~vtkGDCMMedicalImageProperties	. 877
27.352.2Member Function Documentation	. 877
27.352.2.1Clear	. 877
27.352.2.2GetFile	. 877
27.352.2.3New	. 877
27.352.2.4PrintSelf	. 877
27.352.2.5PushBackFile	. 877
27.352.2.6vtkTypeRevisionMacro	. 877
27.352.3Friends And Related Function Documentation	. 877
27.352.3.1vtkGDCMImageReader	. 877
27.352.3.2vtkGDCMImageReader2	. 877
27.352.3.3vtkGDCMImageWriter	. 877
27.353.1vtkGDCMPolyDataReader Class Reference	. 877
27.353.1.1Detailed Description	. 879
27.353.1.2Constructor & Destructor Documentation	. 879
27.353.1.2.1vtkGDCMPolyDataReader	. 879
27.353.1.2.2~vtkGDCMPolyDataReader	. 879
27.353.1.3Member Function Documentation	. 879
27.353.1.3.1FillMedicalImageInformation	. 879
27.353.1.3.2New	. 879
27.353.1.3.3PrintSelf	. 879
27.353.1.3.4RequestData	. 879
27.353.1.3.5RequestData_HemodynamicWaveformStorage	. 879
27.353.1.3.6RequestData_RTStructureSetStorage	. 880
27.353.1.3.7RequestInformation	. 880
27.353.1.3.8RequestInformation_HemodynamicWaveformStorage	. 880
27.353.1.3.9RequestInformation_RTStructureSetStorage	. 880
27.353.1.3.10GetObjectMacro	. 880
27.353.1.3.11GetObjectMacro	. 880
27.353.1.3.12GetStringMacro	. 880
27.353.1.3.13SetStringMacro	. 880
27.353.1.3.14TypeRevisionMacro	. 880
27.353.1.4Member Data Documentation	. 880
27.353.1.4.1FileName	. 880
27.353.1.4.2MedicalImageProperties	. 880
27.353.1.4.3RTStructSetProperties	. 880
27.354.1vtkGDCMPolyDataWriter Class Reference	. 880

27.354.1	Detailed Description	882
27.354.2	Constructor & Destructor Documentation	882
27.354.2.1	vtkGDCMPolyDataWriter	882
27.354.2.2	~vtkGDCMPolyDataWriter	882
27.354.3	Member Function Documentation	882
27.354.3.1	InitializeRTStructSet	882
27.354.3.2	New	882
27.354.3.3	PrintSelf	882
27.354.3.4	SetMedicalImageProperties	882
27.354.3.5	SetNumberOfInputPorts	883
27.354.3.6	SetRTStructSetProperties	883
27.354.3.7	vtkTypeRevisionMacro	883
27.354.3.8	WriteData	883
27.354.3.9	WriteRTSTRUCTData	883
27.354.3.10	WriteRTSTRUCTInfo	883
27.354.4	Member Data Documentation	883
27.354.4.1	MedicalImageProperties	883
27.354.4.2	RTStructSetProperties	883
27.355	vtkGDCMTesting Class Reference	883
27.355.1	Detailed Description	884
27.355.2	Member Typedef Documentation	885
27.355.2.1	MD5MetalImagesType	885
27.355.3	Constructor & Destructor Documentation	885
27.355.3.1	vtkGDCMTesting	885
27.355.3.2	~vtkGDCMTesting	885
27.355.4	Member Function Documentation	885
27.355.4.1	GetGDCMDataRoot	885
27.355.4.2	GetMD5MetalImage	885
27.355.4.3	GetMHDMD5FromFile	885
27.355.4.4	GetNumberOfMD5MetalImages	885
27.355.4.5	GetRAWMD5FromFile	885
27.355.4.6	GetVTKDataRoot	885
27.355.4.7	New	885
27.355.4.8	PrintSelf	886
27.355.4.9	vtkTypeRevisionMacro	886
27.356	vtkGDCMThreadedImageReader Class Reference	886
27.356.1	Constructor & Destructor Documentation	887

27.356.1.1	vtkGDCMThreadedImageReader	887
27.356.1.2	~vtkGDCMThreadedImageReader	888
27.356.2	Member Function Documentation	888
27.356.2.1	ExecuteData	888
27.356.2.2	ExecuteInformation	888
27.356.2.3	New	888
27.356.2.4	PrintSelf	888
27.356.2.5	ReadFiles	888
27.356.2.6	RequestDataCompat	888
27.356.2.7	vtkBooleanMacro	888
27.356.2.8	vtkGetMacro	888
27.356.2.9	vtkSetMacro	888
27.356.2.10	vtkSetMacro	888
27.356.2.11	vtkSetMacro	888
27.356.2.12	vtkTypeRevisionMacro	888
27.357	vtkGDCMThreadedImageReader2 Class Reference	888
27.357.1	Constructor & Destructor Documentation	890
27.357.1.1	vtkGDCMThreadedImageReader2	890
27.357.1.2	~vtkGDCMThreadedImageReader2	890
27.357.2	Member Function Documentation	890
27.357.2.1	GetFileName	890
27.357.2.2	New	890
27.357.2.3	PrintSelf	890
27.357.2.4	RequestInformation	890
27.357.2.5	SetFileName	890
27.357.2.6	SetFileNames	890
27.357.2.7	SplitExtent	891
27.357.2.8	ThreadedRequestData	891
27.357.2.9	vtkBooleanMacro	891
27.357.2.10	vtkBooleanMacro	891
27.357.2.11	vtkBooleanMacro	891
27.357.2.12	vtkGetMacro	891
27.357.2.13	vtkGetMacro	891
27.357.2.14	vtkGetMacro	891
27.357.2.15	vtkGetMacro	891
27.357.2.16	vtkGetMacro	891
27.357.2.17	vtkGetMacro	891

27.357.2.18	GetMacro	891
27.357.2.19	GetMacro	891
27.357.2.20	GetObjectMacro	891
27.357.2.21	GetVector3Macro	891
27.357.2.22	GetVector3Macro	891
27.357.2.23	GetVector6Macro	891
27.357.2.24	SetMacro	891
27.357.2.25	SetMacro	891
27.357.2.26	SetMacro	891
27.357.2.27	SetMacro	891
27.357.2.28	SetMacro	891
27.357.2.29	SetMacro	891
27.357.2.30	SetMacro	891
27.357.2.31	SetVector3Macro	891
27.357.2.32	SetVector3Macro	892
27.357.2.33	SetVector6Macro	892
27.357.2.34	TypeRevisionMacro	892
27.358	vtkImageColorViewer Class Reference	892
27.358.1	Detailed Description	895
27.358.2	Member Enumeration Documentation	895
27.358.2.1	anonymous enum	895
27.358.3	Constructor & Destructor Documentation	895
27.358.3.1	vtkImageColorViewer	895
27.358.3.2	~vtkImageColorViewer	895
27.358.4	Member Function Documentation	895
27.358.4.1	AddInput	895
27.358.4.2	AddInputConnection	895
27.358.4.3	GetColorLevel	895
27.358.4.4	GetColorWindow	895
27.358.4.5	GetInput	895
27.358.4.6	GetOffScreenRendering	896
27.358.4.7	GetOverlayVisibility	896
27.358.4.8	GetPosition	896
27.358.4.9	GetSize	896
27.358.4.10	GetSliceMax	896
27.358.4.11	GetSliceMin	896
27.358.4.12	GetSliceRange	896

27.358.4.10	GetSliceRange	. . . . .	896
27.358.4.10	GetSliceRange	. . . . .	896
27.358.4.10	GetWindowName	. . . . .	896
27.358.4.10	InstallPipeline	. . . . .	896
27.358.4.11	New	. . . . .	896
27.358.4.11	PrintSelf	. . . . .	896
27.358.4.11	Render	. . . . .	896
27.358.4.20	SetColorLevel	. . . . .	896
27.358.4.20	SetColorWindow	. . . . .	896
27.358.4.20	SetDisplayId	. . . . .	896
27.358.4.20	SetInput	. . . . .	896
27.358.4.20	SetInputConnection	. . . . .	897
27.358.4.25	SetOffScreenRendering	. . . . .	897
27.358.4.26	SetOverlayVisibility	. . . . .	897
27.358.4.27	SetParentId	. . . . .	897
27.358.4.28	SetPosition	. . . . .	897
27.358.4.29	SetPosition	. . . . .	897
27.358.4.30	SetRenderer	. . . . .	897
27.358.4.30	SetRenderWindow	. . . . .	897
27.358.4.30	SetSize	. . . . .	897
27.358.4.30	SetSize	. . . . .	897
27.358.4.30	SetSlice	. . . . .	897
27.358.4.35	SetSliceOrientation	. . . . .	897
27.358.4.36	SetSliceOrientationToXY	. . . . .	897
27.358.4.37	SetSliceOrientationToXZ	. . . . .	897
27.358.4.38	SetSliceOrientationToYZ	. . . . .	897
27.358.4.39	SetupInteractor	. . . . .	898
27.358.4.40	SetWindowId	. . . . .	898
27.358.4.41	UnInstallPipeline	. . . . .	898
27.358.4.42	UpdateDisplayExtent	. . . . .	898
27.358.4.43	UpdateOrientation	. . . . .	898
27.358.4.44	WTK_LEGACY	. . . . .	898
27.358.4.45	WTK_LEGACY	. . . . .	898
27.358.4.46	WTK_LEGACY	. . . . .	898
27.358.4.47	WTK_LEGACY	. . . . .	898
27.358.4.48	BooleanMacro	. . . . .	898
27.358.4.49	GetMacro	. . . . .	898

27.358.4.50	GetMacro	898
27.358.4.51	GetObjectMacro	898
27.358.4.52	GetObjectMacro	898
27.358.4.53	GetObjectMacro	898
27.358.4.54	GetObjectMacro	898
27.358.4.55	GetObjectMacro	898
27.358.4.56	TypeRevisionMacro	898
27.358.5	Friends And Related Function Documentation	898
27.358.5.1	vtkImageColorViewerCallback	898
27.358.6	Member Data Documentation	898
27.358.6.1	FirstRender	898
27.358.6.2	ImageActor	898
27.358.6.3	Interactor	898
27.358.6.4	InteractorStyle	899
27.358.6.5	OverlayImageActor	899
27.358.6.6	Renderer	899
27.358.6.7	RenderWindow	899
27.358.6.8	Slice	899
27.358.6.9	SliceOrientation	899
27.358.6.10	WindowLevel	899
27.359	vtkImageMapToColors16 Class Reference	899
27.359.1	Constructor & Destructor Documentation	901
27.359.1.1	vtkImageMapToColors16	901
27.359.1.2	~vtkImageMapToColors16	901
27.359.2	Member Function Documentation	901
27.359.2.1	GetMTime	901
27.359.2.2	New	901
27.359.2.3	PrintSelf	901
27.359.2.4	RequestData	901
27.359.2.5	RequestInformation	901
27.359.2.6	SetLookupTable	901
27.359.2.7	SetOutputFormatToLuminance	901
27.359.2.8	SetOutputFormatToLuminanceAlpha	901
27.359.2.9	SetOutputFormatToRGB	901
27.359.2.10	SetOutputFormatToRGBA	901
27.359.2.11	ThreadedRequestData	901
27.359.2.12	vtkBooleanMacro	901

27.359.2.1	tkGetMacro	. . . . .	901
27.359.2.1	tkGetMacro	. . . . .	901
27.359.2.1	tkGetMacro	. . . . .	901
27.359.2.1	tkGetObjectMacro	. . . . .	901
27.359.2.1	tkSetMacro	. . . . .	902
27.359.2.1	tkSetMacro	. . . . .	902
27.359.2.1	tkSetMacro	. . . . .	902
27.359.2.2	tkTypeRevisionMacro	. . . . .	902
27.359.3	Member Data Documentation	. . . . .	902
27.359.3.1	ActiveComponent	. . . . .	902
27.359.3.2	DataWasPassed	. . . . .	902
27.359.3.3	LookupTable	. . . . .	902
27.359.3.4	OutputFormat	. . . . .	902
27.359.3.5	PassAlphaToOutput	. . . . .	902
27.360	tkImageMapToWindowLevelColors2 Class Reference	. . . . .	902
27.360	Constructor & Destructor Documentation	. . . . .	904
27.360.1	tkImageMapToWindowLevelColors2	. . . . .	904
27.360.1.2	~tkImageMapToWindowLevelColors2	. . . . .	904
27.360.2	Member Function Documentation	. . . . .	904
27.360.2.1	New	. . . . .	904
27.360.2.2	PrintSelf	. . . . .	904
27.360.2.3	RequestData	. . . . .	904
27.360.2.4	RequestInformation	. . . . .	904
27.360.2.5	ThreadedRequestData	. . . . .	904
27.360.2.6	tkGetMacro	. . . . .	904
27.360.2.7	tkGetMacro	. . . . .	904
27.360.2.8	tkSetMacro	. . . . .	904
27.360.2.9	tkSetMacro	. . . . .	904
27.360.2.10	tkTypeRevisionMacro	. . . . .	904
27.360.3	Member Data Documentation	. . . . .	904
27.360.3.1	Level	. . . . .	904
27.360.3.2	Window	. . . . .	904
27.361	tkImagePlanarComponentsToComponents Class Reference	. . . . .	904
27.361	Constructor & Destructor Documentation	. . . . .	906
27.361.1	tkImagePlanarComponentsToComponents	. . . . .	906
27.361.1.2	~tkImagePlanarComponentsToComponents	. . . . .	906
27.361.2	Member Function Documentation	. . . . .	906



27.361.2.1New	906
27.361.2.2PrintSelf	906
27.361.2.3RequestData	906
27.361.2.4vtkTypeRevisionMacro	906
27.362.vtkImageRGBToYBR Class Reference	906
27.362.1.Constructor & Destructor Documentation	907
27.362.1.1vtkImageRGBToYBR	907
27.362.1.2~vtkImageRGBToYBR	907
27.362.2.Member Function Documentation	907
27.362.2.1New	907
27.362.2.2PrintSelf	907
27.362.2.3ThreadedExecute	907
27.362.2.4vtkTypeRevisionMacro	907
27.363.vtkImageYBRToRGB Class Reference	908
27.363.1.Constructor & Destructor Documentation	909
27.363.1.1vtkImageYBRToRGB	909
27.363.1.2~vtkImageYBRToRGB	909
27.363.2.Member Function Documentation	909
27.363.2.1New	909
27.363.2.2PrintSelf	909
27.363.2.3ThreadedExecute	909
27.363.2.4vtkTypeRevisionMacro	909
27.364.vtkLookupTable16 Class Reference	909
27.364.1.Constructor & Destructor Documentation	910
27.364.1.1vtkLookupTable16	910
27.364.1.2~vtkLookupTable16	910
27.364.2.Member Function Documentation	910
27.364.2.1Build	911
27.364.2.2GetPointer	911
27.364.2.3MapScalarsThroughTable2	911
27.364.2.4New	911
27.364.2.5PrintSelf	911
27.364.2.6SetNumberOfTableValues	911
27.364.2.7vtkTypeRevisionMacro	911
27.364.2.8WritePointer	911
27.364.3.Member Data Documentation	911
27.364.3.1Table16	911

27.365.1	vtkRTStructSetProperties Class Reference	911
27.365.1	Detailed Description	913
27.365.2	Constructor & Destructor Documentation	913
27.365.2.1	vtkRTStructSetProperties	913
27.365.2.2	~vtkRTStructSetProperties	913
27.365.3	Member Function Documentation	913
27.365.3.1	AddContourReferencedFrameOfReference	914
27.365.3.2	AddReferencedFrameOfReference	914
27.365.3.3	AddStructureSetROI	914
27.365.3.4	AddStructureSetROIObservation	914
27.365.3.5	Clear	914
27.365.3.6	DeepCopy	914
27.365.3.7	GetContourReferencedFrameOfReferenceClassUID	914
27.365.3.8	GetContourReferencedFrameOfReferenceInstanceUID	914
27.365.3.9	GetNumberOfContourReferencedFrameOfReferences	914
27.365.3.10	GetNumberOfContourReferencedFrameOfReferences	914
27.365.3.11	GetNumberOfReferencedFrameOfReferences	914
27.365.3.12	GetNumberOfStructureSetROIs	914
27.365.3.13	GetReferencedFrameOfReferenceClassUID	914
27.365.3.14	GetReferencedFrameOfReferenceInstanceUID	914
27.365.3.15	SetStructureSetObservationNumber	914
27.365.3.16	GetStructureSetROIDescription	914
27.365.3.17	GetStructureSetROIGenerationAlgorithm	914
27.365.3.18	GetStructureSetROIName	914
27.365.3.19	GetStructureSetROINumber	914
27.365.3.20	GetStructureSetROIObservationLabel	914
27.365.3.21	GetStructureSetROIRefFrameRefUID	914
27.365.3.22	GetStructureSetRTROIInterpretedType	914
27.365.3.23	New	915
27.365.3.24	PrintSelf	915
27.365.3.25	GetStringMacro	915
27.365.3.26	GetStringMacro	915
27.365.3.27	GetStringMacro	915
27.365.3.28	GetStringMacro	915
27.365.3.29	GetStringMacro	915
27.365.3.30	GetStringMacro	915
27.365.3.31	GetStringMacro	915

27.365.3.32kGetStringMacro . . . . .	915
27.365.3.33kGetStringMacro . . . . .	915
27.365.3.34kSetStringMacro . . . . .	915
27.365.3.35kSetStringMacro . . . . .	915
27.365.3.36kSetStringMacro . . . . .	915
27.365.3.37kSetStringMacro . . . . .	915
27.365.3.38kSetStringMacro . . . . .	915
27.365.3.39kSetStringMacro . . . . .	915
27.365.3.40kSetStringMacro . . . . .	915
27.365.3.41kSetStringMacro . . . . .	915
27.365.3.42kSetStringMacro . . . . .	915
27.365.3.43kTypeRevisionMacro . . . . .	915
27.365.4Member Data Documentation . . . . .	915
27.365.4.1Internals . . . . .	915
27.365.4.2ReferenceFrameOfReferenceUID . . . . .	915
27.365.4.3ReferenceSeriesInstanceUID . . . . .	916
27.365.4.4SeriesInstanceUID . . . . .	916
27.365.4.5SOPInstanceUID . . . . .	916
27.365.4.6StructureSetDate . . . . .	916
27.365.4.7StructureSetLabel . . . . .	916
27.365.4.8StructureSetName . . . . .	916
27.365.4.9StructureSetTime . . . . .	916
27.365.4.10StudyInstanceUID . . . . .	916
27.366dcm::Waveform Class Reference . . . . .	916
27.366.1Detailed Description . . . . .	916
27.366.2Constructor & Destructor Documentation . . . . .	916
27.366.2.1Waveform . . . . .	916
27.367dcm::Writer Class Reference . . . . .	916
27.367.1Detailed Description . . . . .	919
27.367.2Constructor & Destructor Documentation . . . . .	920
27.367.2.1Writer . . . . .	920
27.367.2.2~Writer . . . . .	920
27.367.3Member Function Documentation . . . . .	920
27.367.3.1CheckFileMetaInformationOff . . . . .	920
27.367.3.2CheckFileMetaInformationOn . . . . .	920
27.367.3.3GetFile . . . . .	920
27.367.3.4GetStreamPtr . . . . .	920

27.367.3.5SetCheckFileMetaInformation . . . . .	920
27.367.3.6SetFile . . . . .	920
27.367.3.7SetFileName . . . . .	920
27.367.3.8SetStream . . . . .	921
27.367.3.9SetWriteDataSetOnly . . . . .	921
27.367.3.10Write . . . . .	921
27.367.4Friends And Related Function Documentation . . . . .	921
27.367.4.1StreamImageWriter . . . . .	921
27.367.5Member Data Documentation . . . . .	921
27.367.5.1Ostream . . . . .	921
27.367.5.2SStream . . . . .	921
27.368dcm::XMLDictReader Class Reference . . . . .	921
27.368.1Detailed Description . . . . .	922
27.368.2Constructor & Destructor Documentation . . . . .	923
27.368.2.1XMLDictReader . . . . .	923
27.368.2.2~XMLDictReader . . . . .	923
27.368.3Member Function Documentation . . . . .	923
27.368.3.1CharacterDataHandler . . . . .	923
27.368.3.2EndElement . . . . .	923
27.368.3.3GetDict . . . . .	923
27.368.3.4HandleDescription . . . . .	923
27.368.3.5HandleEntry . . . . .	923
27.368.3.6StartElement . . . . .	923
27.369dcm::XMLPrinter Class Reference . . . . .	923
27.369.1Member Enumeration Documentation . . . . .	925
27.369.1.1PrintStyles . . . . .	925
27.369.2Constructor & Destructor Documentation . . . . .	925
27.369.2.1XMLPrinter . . . . .	925
27.369.2.2~XMLPrinter . . . . .	925
27.369.3Member Function Documentation . . . . .	925
27.369.3.1GetPrintStyle . . . . .	925
27.369.3.2HandleBulkData . . . . .	925
27.369.3.3Print . . . . .	925
27.369.3.4PrintDataElement . . . . .	925
27.369.3.5PrintDataSet . . . . .	925
27.369.3.6PrintSQ . . . . .	925
27.369.3.7SetFile . . . . .	925

27.369.3.8SetStyle . . . . .	925
27.369.4Member Data Documentation . . . . .	925
27.369.4.1F . . . . .	925
27.369.4.2PrintStyle . . . . .	925
27.370gdcm::XMLPrivateDictReader Class Reference . . . . .	926
27.370.1Detailed Description . . . . .	927
27.370.2Constructor & Destructor Documentation . . . . .	927
27.370.2.1XMLPrivateDictReader . . . . .	927
27.370.2.2~XMLPrivateDictReader . . . . .	927
27.370.3Member Function Documentation . . . . .	927
27.370.3.1CharacterDataHandler . . . . .	927
27.370.3.2EndElement . . . . .	927
27.370.3.3GetPrivateDict . . . . .	927
27.370.3.4HandleDescription . . . . .	927
27.370.3.5HandleEntry . . . . .	927
27.370.3.6StartElement . . . . .	927
<b>28 File Documentation</b>	<b>929</b>
28.1 gdcm2pnm.dox File Reference . . . . .	929
28.2 gdcm2vtk.dox File Reference . . . . .	929
28.3 gdcmAAbortPDU.h File Reference . . . . .	929
28.4 gdcmAAssociateACPDU.h File Reference . . . . .	930
28.5 gdcmAAssociateRJPDU.h File Reference . . . . .	930
28.6 gdcmAAssociateRQPDU.h File Reference . . . . .	931
28.7 gdcmAbstractSyntax.h File Reference . . . . .	932
28.8 gdcmanon.dox File Reference . . . . .	933
28.9 gdcmAnonymizeEvent.h File Reference . . . . .	933
28.10gdcmAnonymizer.h File Reference . . . . .	934
28.11gdcmApplicationContext.h File Reference . . . . .	935
28.12gdcmApplicationEntity.h File Reference . . . . .	936
28.13gdcmAReleaseRPPDU.h File Reference . . . . .	936
28.14gdcmAReleaseRQPDU.h File Reference . . . . .	937
28.15gdcmARTIMTimer.h File Reference . . . . .	938
28.16gdcmASN1.h File Reference . . . . .	939
28.17gdcmAsynchronousOperationsWindowSub.h File Reference . . . . .	939
28.18gdcmAttribute.h File Reference . . . . .	940
28.19gdcmAudioCodec.h File Reference . . . . .	942

28.20gdcmbase64.h File Reference . . . . .	942
28.21gdcmbaseCompositeMessage.h File Reference . . . . .	943
28.22gdcmbasePDU.h File Reference . . . . .	944
28.23gdcmbaseRootQuery.h File Reference . . . . .	945
28.24gdcmbasicOffsetTable.h File Reference . . . . .	946
28.25gdcmbitmap.h File Reference . . . . .	948
28.26gdcmbitmapToBitmapFilter.h File Reference . . . . .	949
28.27gdcmboxRegion.h File Reference . . . . .	949
28.28gdcmbyteBuffer.h File Reference . . . . .	950
28.29gdcmbyteSwap.h File Reference . . . . .	951
28.30gdcmbyteSwapFilter.h File Reference . . . . .	952
28.31gdcmbyteValue.h File Reference . . . . .	953
28.32gdcmcapiCryptoFactory.h File Reference . . . . .	954
28.33gdcmcapiCryptographicMessageSyntax.h File Reference . . . . .	954
28.34gdcmcEchoMessages.h File Reference . . . . .	955
28.35gdcmcFindMessages.h File Reference . . . . .	956
28.36gdcmcMoveMessages.h File Reference . . . . .	957
28.37gdcmcCodec.h File Reference . . . . .	957
28.38gdcmcCoder.h File Reference . . . . .	958
28.39gdcmcCodeString.h File Reference . . . . .	959
28.40gdcmcCommand.h File Reference . . . . .	960
28.41gdcmcCommandDataSet.h File Reference . . . . .	962
28.42gdcmcCompositeMessageFactory.h File Reference . . . . .	962
28.43gdcmcCompositeNetworkFunctions.h File Reference . . . . .	963
28.44gdcmcConstCharWrapper.h File Reference . . . . .	964
28.45gdcmcconv.dox File Reference . . . . .	964
28.46gdcmcCP246ExplicitDataElement.h File Reference . . . . .	965
28.47gdcmcCryptoFactory.h File Reference . . . . .	965
28.48gdcmcCryptographicMessageSyntax.h File Reference . . . . .	966
28.49gdcmcCSAElement.h File Reference . . . . .	967
28.50gdcmcCSAHeader.h File Reference . . . . .	968
28.51gdcmcCSAHeaderDict.h File Reference . . . . .	969
28.52gdcmcCSAHeaderDictEntry.h File Reference . . . . .	971
28.53gdcmcCStoreMessages.h File Reference . . . . .	972
28.54gdcmcCurve.h File Reference . . . . .	973
28.55gdcmcDataElement.h File Reference . . . . .	974
28.56gdcmcDataEvent.h File Reference . . . . .	975

28.57gdcmlDataSet.h File Reference . . . . .	976
28.58gdcmlDataSetEvent.h File Reference . . . . .	977
28.59gdcmlDataSetHelper.h File Reference . . . . .	978
28.60gdcmlDecoder.h File Reference . . . . .	979
28.61gdcmlDefinedTerms.h File Reference . . . . .	980
28.62gdcmlDeflateStream.h File Reference . . . . .	981
28.63gdcmlDefs.h File Reference . . . . .	981
28.64gdcmlDeltaEncodingCodec.h File Reference . . . . .	983
28.65gdcmlDICOMDIR.h File Reference . . . . .	983
28.66gdcmlDICOMDIRGenerator.h File Reference . . . . .	984
28.67gdcmlDict.h File Reference . . . . .	985
28.68gdcmlDictConverter.h File Reference . . . . .	987
28.69gdcmlDictEntry.h File Reference . . . . .	987
28.70gdcmlDictPrinter.h File Reference . . . . .	989
28.71gdcmlDicts.h File Reference . . . . .	989
28.72gdcmlDiff.dox File Reference . . . . .	990
28.73gdcmlDIMSE.h File Reference . . . . .	991
28.74gdcmlDirectionCosines.h File Reference . . . . .	991
28.75gdcmlDirectory.h File Reference . . . . .	992
28.76gdcmlDirectoryHelper.h File Reference . . . . .	993
28.77gdcmlDummyValueGenerator.h File Reference . . . . .	994
28.78gdcmlDump.dox File Reference . . . . .	994
28.79gdcmlDumper.h File Reference . . . . .	994
28.80gdcmlElement.h File Reference . . . . .	995
28.81gdcmlEncapsulatedDocument.h File Reference . . . . .	997
28.82gdcmlEnumeratedValues.h File Reference . . . . .	997
28.83gdcmlEvent.h File Reference . . . . .	998
28.83.1 Macro Definition Documentation . . . . .	999
28.83.1.1 gdcmlEventMacro . . . . .	999
28.84gdcmlException.h File Reference . . . . .	1000
28.85gdcmlExplicitDataElement.h File Reference . . . . .	1000
28.86gdcmlExplicitImplicitDataElement.h File Reference . . . . .	1001
28.87gdcmlFiducials.h File Reference . . . . .	1002
28.88gdcmlFile.h File Reference . . . . .	1003
28.89gdcmlFileAnonymizer.h File Reference . . . . .	1004
28.90gdcmlFileChangeTransferSyntax.h File Reference . . . . .	1004
28.91gdcmlFileDerivation.h File Reference . . . . .	1005

28.92gdcmlFileExplicitFilter.h File Reference . . . . .	1006
28.93gdcmlFileMetaInformation.h File Reference . . . . .	1007
28.94gdcmlFilename.h File Reference . . . . .	1008
28.95gdcmlFileNameEvent.h File Reference . . . . .	1009
28.96gdcmlFilenameGenerator.h File Reference . . . . .	1010
28.97gdcmlFileSet.h File Reference . . . . .	1010
28.98gdcmlFileStreamer.h File Reference . . . . .	1012
28.99gdcmlFindPatientRootQuery.h File Reference . . . . .	1012
28.100gdcmlFindStudyRootQuery.h File Reference . . . . .	1013
28.101gdcmlFragment.h File Reference . . . . .	1014
28.102gdcmlgendir.dox File Reference . . . . .	1016
28.103gdcmlGlobal.h File Reference . . . . .	1016
28.104gdcmlGroupDict.h File Reference . . . . .	1017
28.105gdcmlIconImage.h File Reference . . . . .	1017
28.106gdcmlIconImageFilter.h File Reference . . . . .	1018
28.107gdcmlIconImageGenerator.h File Reference . . . . .	1019
28.108gdcmlImage.h File Reference . . . . .	1020
28.109gdcmlImageApplyLookupTable.h File Reference . . . . .	1021
28.110gdcmlImageChangePhotometricInterpretation.h File Reference . . . . .	1021
28.111gdcmlImageChangePlanarConfiguration.h File Reference . . . . .	1022
28.112gdcmlImageChangeTransferSyntax.h File Reference . . . . .	1023
28.113gdcmlImageCodec.h File Reference . . . . .	1024
28.114gdcmlImageConverter.h File Reference . . . . .	1025
28.115gdcmlImageFragmentSplitter.h File Reference . . . . .	1026
28.116gdcmlImageHelper.h File Reference . . . . .	1027
28.117gdcmlImageReader.h File Reference . . . . .	1028
28.118gdcmlImageRegionReader.h File Reference . . . . .	1028
28.119gdcmlImageToImageFilter.h File Reference . . . . .	1029
28.120gdcmlImageWriter.h File Reference . . . . .	1030
28.121gdcmlimg.dox File Reference . . . . .	1031
28.122gdcmlImplementationClassUIDSub.h File Reference . . . . .	1031
28.123gdcmlImplementationUIDSub.h File Reference . . . . .	1032
28.124gdcmlImplementationVersionNameSub.h File Reference . . . . .	1033
28.125gdcmlImplicitDataElement.h File Reference . . . . .	1034
28.126gdcmlinfo.dox File Reference . . . . .	1034
28.127gdcmlIOD.h File Reference . . . . .	1035
28.128gdcmlIODEntry.h File Reference . . . . .	1036



28.129	dcmIODs.h File Reference . . . . .	1038
28.130	dcmIPPSorter.h File Reference . . . . .	1039
28.131	dcmItem.h File Reference . . . . .	1040
28.132	dcmJPEG12Codec.h File Reference . . . . .	1042
28.133	dcmJPEG16Codec.h File Reference . . . . .	1042
28.134	dcmJPEG2000Codec.h File Reference . . . . .	1043
28.135	dcmJPEG8Codec.h File Reference . . . . .	1044
28.136	dcmJPEGCodec.h File Reference . . . . .	1045
28.137	dcmJPEGLSCodec.h File Reference . . . . .	1046
28.138	dcmJSON.h File Reference . . . . .	1047
28.139	dcmKAKADUCodec.h File Reference . . . . .	1048
28.140	dcmLegacyMacro.h File Reference . . . . .	1048
28.140.1	Macro Definition Documentation . . . . .	1049
28.140.1.1	1GDCM_LEGACY . . . . .	1049
28.140.1.2	2GDCM_LEGACY_BODY . . . . .	1049
28.140.1.3	3GDCM_LEGACY_REPLACED_BODY . . . . .	1049
28.141	dcmLO.h File Reference . . . . .	1049
28.142	dcmLookupTable.h File Reference . . . . .	1050
28.143	dcmMacro.h File Reference . . . . .	1051
28.144	dcmMacroEntry.h File Reference . . . . .	1053
28.144.1	Macro Definition Documentation . . . . .	1054
28.144.1.1	1GDCMMACROENTRY_H . . . . .	1054
28.145	dcmMacros.h File Reference . . . . .	1054
28.146	dcmMaximumLengthSub.h File Reference . . . . .	1056
28.147	dcmMD5.h File Reference . . . . .	1057
28.148	dcmMediaStorage.h File Reference . . . . .	1058
28.149	dcmMeshPrimitive.h File Reference . . . . .	1059
28.150	dcmModule.h File Reference . . . . .	1060
28.151	dcmModuleEntry.h File Reference . . . . .	1062
28.152	dcmModules.h File Reference . . . . .	1063
28.153	dcmMovePatientRootQuery.h File Reference . . . . .	1065
28.154	dcmMoveStudyRootQuery.h File Reference . . . . .	1066
28.155	dcmNestedModuleEntries.h File Reference . . . . .	1066
28.156	dcmNetworkEvents.h File Reference . . . . .	1068
28.157	dcmNetworkStateID.h File Reference . . . . .	1069
28.158	dcmObject.h File Reference . . . . .	1070
28.159	dcmOpenSSLCryptoFactory.h File Reference . . . . .	1070

28.160dcmOpenSSLCryptographicMessageSyntax.h File Reference . . . . .	1071
28.164dcmOpenSSLP7CryptoFactory.h File Reference . . . . .	1072
28.162dcmOpenSSLP7CryptographicMessageSyntax.h File Reference . . . . .	1073
28.163dcmOrientation.h File Reference . . . . .	1075
28.164dcmOverlay.h File Reference . . . . .	1075
28.165dcmpap3.dox File Reference . . . . .	1076
28.166dcmParseException.h File Reference . . . . .	1076
28.167dcmParser.h File Reference . . . . .	1078
28.168dcmPatient.h File Reference . . . . .	1078
28.169dcmPDDataTFPDU.h File Reference . . . . .	1079
28.170dcmPDBElement.h File Reference . . . . .	1080
28.174dcmPDBHeader.h File Reference . . . . .	1082
28.172dcmpdf.dox File Reference . . . . .	1082
28.173dcmPDFCodec.h File Reference . . . . .	1082
28.174dcmPDUFactory.h File Reference . . . . .	1083
28.175dcmPersonName.h File Reference . . . . .	1084
28.176dcmPGXCodec.h File Reference . . . . .	1084
28.177dcmPhotometricInterpretation.h File Reference . . . . .	1085
28.178dcmPixelFormat.h File Reference . . . . .	1086
28.179dcmPixmap.h File Reference . . . . .	1087
28.180dcmPixmapReader.h File Reference . . . . .	1088
28.184dcmPixmapToPixmapFilter.h File Reference . . . . .	1089
28.182dcmPixmapWriter.h File Reference . . . . .	1090
28.183dcmPNMCodec.h File Reference . . . . .	1091
28.184dcmPreamble.h File Reference . . . . .	1092
28.185dcmPresentationContext.h File Reference . . . . .	1093
28.186dcmPresentationContextAC.h File Reference . . . . .	1094
28.187dcmPresentationContextGenerator.h File Reference . . . . .	1095
28.189dcmPresentationContextRQ.h File Reference . . . . .	1096
28.189dcmPresentationDataValue.h File Reference . . . . .	1097
28.190dcmPrinter.h File Reference . . . . .	1097
28.194dcmPrivateTag.h File Reference . . . . .	1098
28.192dcmProgressEvent.h File Reference . . . . .	1100
28.193dcmPVRGCodec.h File Reference . . . . .	1100
28.194dcmPythonFilter.h File Reference . . . . .	1101
28.195dcmQueryBase.h File Reference . . . . .	1102
28.196dcmQueryFactory.h File Reference . . . . .	1103

28.197dcmQueryImage.h File Reference . . . . .	1104
28.198dcmQueryPatient.h File Reference . . . . .	1105
28.199dcmQuerySeries.h File Reference . . . . .	1106
28.200dcmQueryStudy.h File Reference . . . . .	1107
28.201dcmraw.dox File Reference . . . . .	1107
28.202dcmRAWCodec.h File Reference . . . . .	1107
28.203dcmReader.h File Reference . . . . .	1108
28.204dcmRegion.h File Reference . . . . .	1110
28.205dcmRescaler.h File Reference . . . . .	1111
28.206dcmRLECodec.h File Reference . . . . .	1111
28.207dcmRoleSelectionSub.h File Reference . . . . .	1112
28.208dcmscanner.dox File Reference . . . . .	1113
28.209dcmScanner.h File Reference . . . . .	1113
28.210dcmscu.dox File Reference . . . . .	1114
28.211dcmSegment.h File Reference . . . . .	1114
28.212dcmSegmentedPaletteColorLookupTable.h File Reference . . . . .	1115
28.213dcmSegmentHelper.h File Reference . . . . .	1116
28.214dcmSegmentReader.h File Reference . . . . .	1117
28.215dcmSegmentWriter.h File Reference . . . . .	1118
28.216dcmSequenceOfFragments.h File Reference . . . . .	1119
28.217dcmSequenceOfItems.h File Reference . . . . .	1120
28.218dcmSerieHelper.h File Reference . . . . .	1121
28.219dcmSeries.h File Reference . . . . .	1122
28.220dcmServiceClassApplicationInformation.h File Reference . . . . .	1124
28.221dcmServiceClassUser.h File Reference . . . . .	1125
28.222dcmSHA1.h File Reference . . . . .	1125
28.223dcmSimpleSubjectWatcher.h File Reference . . . . .	1126
28.224dcmSmartPointer.h File Reference . . . . .	1127
28.225dcmSOPClassExtendedNegociationSub.h File Reference . . . . .	1128
28.226dcmSOPClassUIDToIOD.h File Reference . . . . .	1129
28.227dcmSorter.h File Reference . . . . .	1130
28.228dcmSpacing.h File Reference . . . . .	1132
28.229dcmSpectroscopy.h File Reference . . . . .	1132
28.230dcmSplitMosaicFilter.h File Reference . . . . .	1133
28.231dcmStaticAssert.h File Reference . . . . .	1134
28.231.1Macro Definition Documentation . . . . .	1134
28.231.1.1GDCM_DO_JOIN . . . . .	1134

28.231.1.2GDCM_DO_JOIN2 . . . . .	1134
28.231.1.3GDCM_JOIN . . . . .	1134
28.231.1.4GDCM_STATIC_ASSERT . . . . .	1135
28.232dcmStreamImageReader.h File Reference . . . . .	1135
28.233dcmStreamImageWriter.h File Reference . . . . .	1135
28.234dcmString.h File Reference . . . . .	1136
28.235dcmStringFilter.h File Reference . . . . .	1137
28.236dcmStudy.h File Reference . . . . .	1138
28.237dcmSubject.h File Reference . . . . .	1139
28.238dcmSurface.h File Reference . . . . .	1139
28.239dcmSurfaceHelper.h File Reference . . . . .	1141
28.240dcmSurfaceReader.h File Reference . . . . .	1141
28.241dcmSurfaceWriter.h File Reference . . . . .	1142
28.242dcmSwapCode.h File Reference . . . . .	1143
28.243dcmSwapper.h File Reference . . . . .	1144
28.244dcmSystem.h File Reference . . . . .	1145
28.245dcmTable.h File Reference . . . . .	1146
28.246dcmTableEntry.h File Reference . . . . .	1146
28.247dcmTableReader.h File Reference . . . . .	1147
28.248dcmTag.h File Reference . . . . .	1149
28.249dcmTagPath.h File Reference . . . . .	1149
28.250dcmTagToVR.h File Reference . . . . .	1150
28.251dcmTar.dox File Reference . . . . .	1151
28.252dcmTerminal.h File Reference . . . . .	1151
28.253dcmTestDriver.h File Reference . . . . .	1153
28.254dcmTesting.h File Reference . . . . .	1153
28.255dcmTrace.h File Reference . . . . .	1154
28.255.1Macro Definition Documentation . . . . .	1155
28.255.1.1GDCM_FUNCTION . . . . .	1155
28.255.1.2dcmAssertAlwaysMacro . . . . .	1155
28.255.1.3dcmAssertMacro . . . . .	1155
28.255.1.4dcmDebugMacro . . . . .	1156
28.255.1.5dcmErrorMacro . . . . .	1156
28.255.1.6dcmWarningMacro . . . . .	1156
28.256dcmTransferSyntax.h File Reference . . . . .	1157
28.257dcmTransferSyntaxSub.h File Reference . . . . .	1158
28.258dcmType.h File Reference . . . . .	1159

28.259dcmTypes.h File Reference . . . . .	1160
28.260dcmUIDGenerator.h File Reference . . . . .	1161
28.261dcmUIDs.h File Reference . . . . .	1161
28.262dcmULAction.h File Reference . . . . .	1162
28.263dcmULActionAA.h File Reference . . . . .	1163
28.264dcmULActionAE.h File Reference . . . . .	1164
28.265dcmULActionAR.h File Reference . . . . .	1165
28.266dcmULActionDT.h File Reference . . . . .	1166
28.267dcmULBasicCallback.h File Reference . . . . .	1166
28.268dcmULConnection.h File Reference . . . . .	1167
28.269dcmULConnectionCallback.h File Reference . . . . .	1168
28.270dcmULConnectionInfo.h File Reference . . . . .	1169
28.271dcmULConnectionManager.h File Reference . . . . .	1171
28.272dcmULEvent.h File Reference . . . . .	1171
28.273dcmULTransitionTable.h File Reference . . . . .	1172
28.274dcmULWritingCallback.h File Reference . . . . .	1174
28.275dcmUNExplicitDataElement.h File Reference . . . . .	1174
28.276dcmUNExplicitImplicitDataElement.h File Reference . . . . .	1175
28.277dcmUnpacker12Bits.h File Reference . . . . .	1175
28.278dcmUsage.h File Reference . . . . .	1176
28.279dcmUserInformation.h File Reference . . . . .	1178
28.280dcmUUIDGenerator.h File Reference . . . . .	1179
28.281dcmValidate.h File Reference . . . . .	1179
28.282dcmValue.h File Reference . . . . .	1180
28.283dcmValueIO.h File Reference . . . . .	1181
28.284dcmVersion.h File Reference . . . . .	1182
28.285dcmviewer.dox File Reference . . . . .	1183
28.286dcmVL.h File Reference . . . . .	1183
28.287dcmVM.h File Reference . . . . .	1184
28.287.1Macro Definition Documentation . . . . .	1185
28.287.1.1TYPETOLENGTH . . . . .	1185
28.288dcmVR.h File Reference . . . . .	1185
28.288.1Macro Definition Documentation . . . . .	1187
28.288.1.1TYPETOENCODING . . . . .	1187
28.288.1.2VRTypeTemplateCase . . . . .	1187
28.289dcmVR16ExplicitDataElement.h File Reference . . . . .	1187
28.290dcmWaveform.h File Reference . . . . .	1188

28.290	gdcmWin32.h File Reference . . . . .	1188
28.291	Macro Definition Documentation . . . . .	1188
28.291.1	1.GDCM_EXPORT . . . . .	1188
28.292	gdcmWriter.h File Reference . . . . .	1189
28.293	gdcmxml.dox File Reference . . . . .	1190
28.294	gdcmXMLDictReader.h File Reference . . . . .	1190
28.295	gdcmXMLPrinter.h File Reference . . . . .	1190
28.296	gdcmXMLPrivateDictReader.h File Reference . . . . .	1191
28.297	README.txt File Reference . . . . .	1192
28.298	testsList.txt File Reference . . . . .	1192
28.299	tkGDCMImageReader.h File Reference . . . . .	1192
28.299.1	Macro Definition Documentation . . . . .	1193
28.299.1.1	1.VTK_CMYK . . . . .	1193
28.299.1.2	2.VTK_INVERSE_LUMINANCE . . . . .	1193
28.299.1.3	3.VTK_LOOKUP_TABLE . . . . .	1193
28.299.1.4	4.VTK_YBR . . . . .	1193
28.300	tkGDCMImageReader2.h File Reference . . . . .	1193
28.300.1	Macro Definition Documentation . . . . .	1194
28.300.1.1	1.VTK_CMYK . . . . .	1194
28.300.1.2	2.VTK_INVERSE_LUMINANCE . . . . .	1194
28.300.1.3	3.VTK_LOOKUP_TABLE . . . . .	1194
28.300.1.4	4.VTK_YBR . . . . .	1194
28.301	tkGDCMImageWriter.h File Reference . . . . .	1194
28.302	tkGDCMMedicalImageProperties.h File Reference . . . . .	1194
28.303	tkGDCMPolyDataReader.h File Reference . . . . .	1195
28.304	tkGDCMPolyDataWriter.h File Reference . . . . .	1196
28.305	tkGDCMTesting.h File Reference . . . . .	1196
28.306	tkGDCMThreadedImageReader.h File Reference . . . . .	1197
28.307	tkGDCMThreadedImageReader2.h File Reference . . . . .	1198
28.308	tkImageColorViewer.h File Reference . . . . .	1198
28.309	tkImageMapToColors16.h File Reference . . . . .	1199
28.310	tkImageMapToWindowLevelColors2.h File Reference . . . . .	1199
28.311	tkImagePlanarComponentsToComponents.h File Reference . . . . .	1200
28.312	tkImageRGBToYBR.h File Reference . . . . .	1200
28.313	tkImageYBRToRGB.h File Reference . . . . .	1201
28.314	tkLookupTable16.h File Reference . . . . .	1201
28.315	tkRTStructSetProperties.h File Reference . . . . .	1202

<b>29 Example Documentation</b>	<b>1203</b>
29.1 AWTMedical3.java	1203
29.2 BasicAnonymizer.cs	1207
29.3 BasicImageAnonymizer.cs	1208
29.4 CastConvertPhilips.py	1210
29.5 ChangePrivateTags.cxx	1212
29.6 ChangeSequenceUltrasound.cxx	1213
29.7 CheckBigEndianBug.cxx	1214
29.8 ClinicalTrialAnnotate.cxx	1216
29.9 ClinicalTrialIdentificationWorkflow.cs	1217
29.10 CompressImage.cxx	1220
29.11 CompressLossyJPEG.cs	1221
29.12 Compute3DSpacing.cxx	1222
29.13 Convert16BitsTo8Bits.cxx	1223
29.14 ConvertMPL.py	1224
29.15 ConvertMultiFrameToSingleFrame.cxx	1225
29.16 ConvertNumpy.py	1226
29.17 ConvertPIL.py	1227
29.18 ConvertRGBToLuminance.cxx	1229
29.19 ConvertSingleBitTo8Bits.cxx	1229
29.20 ConvertToQImage.cxx	1231
29.21 CreateARGBImage.cxx	1232
29.22 CreateCMYKImage.cxx	1233
29.23 CreateFakeRTDOSE.cxx	1234
29.24 CreateJPIPDataSet.cxx	1236
29.25 CreateRAWStorage.py	1237
29.26 csa2img.cxx	1239
29.27 CStoreQtProgress.cxx	1241
29.28 DecompressImage.cs	1243
29.29 DecompressImage.java	1245
29.30 DecompressImage.py	1246
29.31 DecompressImageMultiframe.cs	1246
29.32 DecompressJPEGFile.cs	1248
29.33 DecompressPixmap.java	1249
29.34 DiffFile.cxx	1250
29.35 DiscriminateVolume.cxx	1251
29.36 DumbAnonymizer.py	1255

29.37DumpADAC.cxx . . . . .	1257
29.38DumpExamCard.cxx . . . . .	1261
29.39DumpGEMSMovieGroup.cxx . . . . .	1268
29.40DumpImageHeaderInfo.cxx . . . . .	1274
29.41DumpPhilipsECHO.cxx . . . . .	1277
29.42DumpToSQLITE3.cxx . . . . .	1282
29.43DuplicatePCDE.cxx . . . . .	1284
29.44ELSCINT1WaveToText.cxx . . . . .	1286
29.45EncapsulateFileInRawData.cxx . . . . .	1288
29.46ExtractEncapsulatedFile.cs . . . . .	1289
29.47ExtractEncryptedContent.cxx . . . . .	1290
29.48ExtractIconFromFile.cxx . . . . .	1291
29.49ExtractImageRegion.cs . . . . .	1292
29.50ExtractImageRegion.java . . . . .	1294
29.51ExtractImageRegionWithLUT.cs . . . . .	1295
29.52Extracting_All_Resolution.cxx . . . . .	1296
29.53ExtractOneFrame.cs . . . . .	1302
29.54Fake_Image_Using_Stream_Image_Writer.cxx . . . . .	1303
29.55FileAnonymize.cs . . . . .	1306
29.56FileAnonymize.java . . . . .	1307
29.57FileChangeTS.cs . . . . .	1308
29.58FileChangeTSLossy.cs . . . . .	1310
29.59FileStreaming.cs . . . . .	1313
29.60FindAllPatientName.py . . . . .	1314
29.61FixBrokenJ2K.cxx . . . . .	1314
29.62FixCommaBug.py . . . . .	1316
29.63FixJAIBugJPEGLS.cxx . . . . .	1317
29.64gdcmmorthoplanes.cxx . . . . .	1320
29.65gdcmmreslice.cxx . . . . .	1326
29.66gdcmmrtionplan.cxx . . . . .	1328
29.67gdcmmrtplan.cxx . . . . .	1333
29.68gdcmmscene.cxx . . . . .	1337
29.69gdcmmtexture.cxx . . . . .	1338
29.70gdcmmvolume.cxx . . . . .	1340
29.71GenAllIVR.cxx . . . . .	1341
29.72GenerateDICOMDIR.cs . . . . .	1344
29.73GenerateRTSTRUCT.cxx . . . . .	1344



29.74GenerateStandardSOPClasses.cxx . . . . .	1347
29.75GenFakeIdentifyFile.cxx . . . . .	1348
29.76GenFakeImage.cxx . . . . .	1351
29.77GenLongSeqs.cxx . . . . .	1352
29.78GenSeqs.cxx . . . . .	1354
29.79GetArray.cs . . . . .	1355
29.80GetJPEGSamplePrecision.cxx . . . . .	1356
29.81GetPortionCSAHeader.py . . . . .	1358
29.82GetSequenceUltrasound.cxx . . . . .	1359
29.83GetSubSequenceData.cxx . . . . .	1361
29.84headsq2dcm.py . . . . .	1363
29.85HelloActiviz.cs . . . . .	1364
29.86HelloActiviz2.cs . . . . .	1365
29.87HelloActiviz3.cs . . . . .	1367
29.88HelloActiviz4.cs . . . . .	1367
29.89HelloActiviz5.cs . . . . .	1368
29.90HelloSimple.java . . . . .	1369
29.91HelloVizWorld.cxx . . . . .	1370
29.92HelloVTKWorld.cs . . . . .	1371
29.93HelloVTKWorld.java . . . . .	1372
29.94HelloVTKWorld2.cs . . . . .	1373
29.95HelloWorld.cxx . . . . .	1374
29.96HelloWorld.py . . . . .	1375
29.97iU22tomultisc.cxx . . . . .	1376
29.98LargeVRDSExplicit.cxx . . . . .	1377
29.99MagnifyFile.cxx . . . . .	1379
29.100ManipulateFile.cs . . . . .	1380
29.101ManipulateFile.py . . . . .	1381
29.102ManipulateSequence.py . . . . .	1383
29.103MergeFile.py . . . . .	1384
29.104MergeTwoFiles.cxx . . . . .	1384
29.105MetalImageMD5Activiz.cs . . . . .	1386
29.106MIPViewer.java . . . . .	1387
29.107MpegVideoInfo.cs . . . . .	1389
29.108MPRViewer.java . . . . .	1394
29.109MPRViewer2.java . . . . .	1396
29.110MrProtocol.cxx . . . . .	1400

29.11NewSequence.cs	1407
29.11NewSequence.py	1408
29.11OffscreenImage.cxx	1409
29.11PatchFile.cxx	1410
29.11PhilipsPrivateRescaleInterceptSlope.py	1412
29.11PlaySound.py	1413
29.11Pmsct_rgb1.cxx	1414
29.11PrivateDict.py	1417
29.11PublicDict.cxx	1418
29.12QIDO-RS.cxx	1419
29.12ReadAndDumpDICOMDIR.cxx	1419
29.12ReadAndDumpDICOMDIR.py	1423
29.12ReadAndPrintAttributes.cxx	1425
29.12ReadExplicitLengthSQIVR.cxx	1427
29.12ReadFiles.java	1427
29.12ReadGEMSSDO.cxx	1428
29.12ReadMultiTimesException.cxx	1431
29.12ReadSeriesIntoVTK.java	1432
29.12ReadUTF8QtDir.cxx	1433
29.13RefCounting.cs	1434
29.13ReformatFile.cs	1435
29.13RemovePrivateTags.py	1436
29.13RescaleImage.cs	1437
29.13Reslicesphere.cxx	1438
29.13ReWriteSCAsMR.py	1446
29.13Re2img.cxx	1447
29.13rstructapp.cxx	1450
29.13ScanDirectory.cs	1451
29.13ScanDirectory.java	1452
29.14ScanDirectory.py	1456
29.14SendFileSCU.cs	1457
29.14SimplePrint.cs	1457
29.14SimplePrintPatientName.cs	1458
29.14SimpleScanner.cxx	1459
29.14SortImage.cxx	1461
29.14SortImage.py	1462
29.14SortImage2.cs	1463

29.148	StandardizeFiles.cs	1464
29.148	StreamImageReaderTest.cxx	1465
29.150	TestByteSwap.cxx	1469
29.151	TestReader.cxx	1471
29.152	TestReader.py	1472
29.153	Threadgdc.cxx	1473
29.154	TraverseModules.cxx	1476
29.155	id_unique.cxx	1477
29.156	VolumeSorter.cxx	1478
29.157	WriteBuffer.py	1480

**Index****1483**



## Chapter 1

# GDCM Documentation

This is the developers documentation.

A PDF version of this doxygen documentation can be found here:

<http://gdcm.sourceforge.net/2.4/gdcm-2.4.4.pdf>

A tarball version of this HTML doxygen documentation can be found here:

<http://gdcm.sourceforge.net/2.4/gdcm-2.4.4-doc.tar.gz>

### Author

Mathieu Malaterre



## Chapter 2

# off-screen rendering of DICOM images

### 2.1 SYNOPSIS

```
gdcm2pnm [options] file-in bitmap-out
```

### 2.2 DESCRIPTION

The **gdcm2pnm** command line program takes as input a DICOM file and produces a rendered bitmap file.

### 2.3 PARAMETERS

file-in    DICOM input filename

bitmap-out    Bitmap output filename

### 2.4 OPTIONS

#### 2.4.1 OPTIONS

#### 2.4.2 general options

```
-h    --help  
      print this help text and exit  
  
-v    --version  
      print version information and exit  
  
-V    --verbose  
      verbose mode (warning+error).  
  
-W    --warning  
      warning mode, print warning information
```

```
-E  --error  
    error mode, print error information  
  
-D  --debug  
    debug mode, print debug information
```

## 2.5 Simple usage

**gdcm2pnm** will take as input DICOM and render it into a bitmap file using the window/level attributes value.

```
$ gdcm2pnm input.dcm output.png
```

It is much different from the **gdcmraw** or **gdcmimg** command line tool as it will render a DICOM image. This means that the output will be rendered in 8bits ready for display.

## 2.6 SEE ALSO

**gdcm2vtk(1)**, **gdcmimg(1)**

## 2.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 3

# Convert a file supported by VTK into DICOM.

### 3.1 SYNOPSIS

```
gdc2vtk [options] file-in file-out
```

### 3.2 DESCRIPTION

The **gdc2vtk** takes as input any file supported by VTK (including DICOM file) and will generate as output a DICOM file.

### 3.3 PARAMETERS

```
file-in    input filename (DICOM or VTK supported)

file-out    output filename (DICOM or VTK supported)
```

### 3.4 OPTIONS

#### 3.4.1 OPTIONS

--force-rescale	force rescale.
--force-spacing	force spacing.
--palette-color	when supported generate a PALETTE COLOR file.
--argb	when supported generate a ARGB file.
--compress	when supported generate a compressed file.
--use-vtkdicom	Use vtkDICOMImageReader (instead of GDCM).
--modality	set Modality.
--lower-left	set lower left.
--shift	set shift.
--scale	set scale.
--compress	set compression (MetaIO).
-T --study-uid	Study UID.
-S --series-uid	Series UID.
--root-uid	Root UID.

### 3.4.2 compression options

```
-J --jpeg          Compress image in jpeg.
-K --j2k          Compress image in j2k.
-L --jpegls       Compress image in jpeg-ls.
-R --rle          Compress image in rle (lossless only).
```

### 3.4.3 general options

```
-h  --help          print this help text and exit
-v  --version       print version information and exit
-V  --verbose       verbose mode (warning+error).
-W  --warning       warning mode, print warning information
-E  --error         error mode, print error information
-D  --debug         debug mode, print debug information
```

### 3.4.4 environment variable

```
GDCM_ROOT_UID Root UID
```

## 3.5 DESCRIPTION

Convert a file supported by VTK into DICOM.

Typical usage is:

```
$ gdcmm2vtk inputfile output.dcm
```

It uses the internal factory mechanism of VTK to recognize a file (CanRead function). See VTK supported file here:

What image file formats can VTK read and write? [http://www.vtk.org/Wiki/VTK\\_FAQ#What\\_image\\_file\\_formats\\_can\\_VTK](http://www.vtk.org/Wiki/VTK_FAQ#What_image_file_formats_can_VTK)

If your input file has 4 components, the 4th comp (alpha) will be removed from the output file as DICOM does not support alpha component anymore (see `--argb` option).

Special care was taken for the following file format:

1. DICOM: Direction Cosines and `vtkMedicalImageInformation` are passed to the output
2. BMP: The file can be saved with a Lookup Table (see `--palette-color`)
3. GE Signa: `vtkMedicalImageProperties` is passed to the output
4. MINC: Direction Cosines is passed to the output
5. TIFF: `vtkTIFFReader` is currently in bad shape in VTK (different behavior in VTK 5.2 and `git/master`). Only u

### 3.5.1 CONVERT MetaImage (mhd, mha)

```
$ gdcmm2vtk inputfile output.mha
```

This command will convert the input DICOM file: inputfile into a MetaImage .mha file. Same goes for .mhd file.

### 3.5.2 CONVERT MHA/MHD

```
$ gdc2vtk inputfile output.mha
```

or

```
$ gdc2vtk inputfile output.mhd
```

This command will convert the input DICOM file: inputfile into a MetaImageData .mha/.mhd file.

### 3.5.3 CONVERT VTI

```
$ gdc2vtk inputfile output.vti
```

This command will convert the input DICOM file: inputfile into a XML VTK ImageData .vti file.

### 3.5.4 CONVERT VTK

```
$ gdc2vtk inputfile output.vtk
```

This command will convert the input DICOM file: inputfile into an old VTK Structured PointSets .vtk file.

## 3.6 CONVERT DICOM

```
$ gdc2vtk input.dcm output.dcm
```

[vtkGDCMImageReader](#) will be used to read in a DICOM file, not the default `vtkDICOMImageReader`. See option `--use-vtkdicom` to use `vtkDICOMImageReader`.

## 3.7 RoundTrip DICOM to MHD to DICOM

```
$ gdc2vtk input_ybr.dcm output.mhd
$ gdc2vtk --modality US --imageformat 7 output.mhd output.dcm
```

The above section shows how to convert a DICOM using the Photometric Interpretation of YBR\_FULL (or even YBR\_FULL\_422 is lossy) into another file format: MetaImage (mhd). Since this file format does not handle color space, we have to explicitly set it using the `--imageformat` command line option. The `--modality` command line option is required in this case since the default Secondary Capture Image Storage Class family does not allow for YBR Photometric Interpretation.

## 3.8 gdc2vtk notes

IMPORTANT NOTE: The internal VTK structured will be filled from the input DICOM, and then pass to the output DICOM writer. Some information might be lost during the conversion DICOM to VTK to DICOM. This option is mostly used to test the `vtkGDCMImageReader/vtkGDCMImageWriter` combination.

IMPORTANT NOTE: When converting from a lossy format such as JPEG, the information of lossiness is important. The output DICOM will contains the required Lossy Image Compression attribute that indicates that image was lossy-compressed somewhere along the pipeline. See also `gdcmimg` (better handling of JPEG in general).

IMPORTANT NOTE: When using `-use-vtkdicom` the output DICOM file will always be written as MR Image Storage as this information is not available from the reader itself. This allow setting the Image Orientation (Patient) properly.

### 3.9 SEE ALSO

`gdcmdump(1)`, `gdcmviewer(1)`, `gdcmimg(1)`

### 3.10 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 4

# Tool to anonymize a DICOM file.

### 4.1 SYNOPSIS

```
gdcmanon [options] file-in file-out
gdcmanon [options] dir-in  dir-out
```

### 4.2 DESCRIPTION

The **gdcmanon** tool is an implementation of PS 3.15 / E.1 / Basic Application Level Confidentiality Profile (Implementation of E.1.1 De-identify & E.1.2 Re-identify)

This tool is split into two very different operating mode:

- An implementation of PS 3.15, see `-e` and `-d` flags
- A dumb mode, see `-dumb`

Dumb mode and PS 3.15 do not work well together, you should really only use one type of anonymization. In case of doubt, avoid using `-dumb`.

In order to use the PS 3.15 implementation (`-d` & `-e` flag), you'll need a certificate to do de-identification operations, and the associated private key to do the re-identification operation. If you are only doing a one-shot anonymization and do not need to properly re-identify the DICOM file, you can safely discard the private key and only keep the certificate. See OpenSSL section below for an example on how to generate the private key/certificate pair.

`gdcmanon` will exit early if OpenSSL was not configured/build properly into the library (see `GDCM_USE_SYSTEM_OPENSSL` in `CMakeLists.txt`).

### 4.3 PARAMETERS

```
file-in  DICOM input filename
```

```
file-out DICOM output filename
```

or

```
file-in  DICOM input directory
```

```
file-out DICOM output directory
```

## 4.4 OPTIONS

You need to specify at least one operating mode, from the following list (and only one):

### 4.4.1 Required parameters

```
-e --de-identify      De-identify DICOM (default)
-d --re-identify      Re-identify DICOM
--dumb                Dumb mode anonymizer
```

Warning when operating in dumb mode, you need to also specify an operation to do, such as 'remove' or 'empty' a tag, see below the dumb mode options.

### 4.4.2 OPTIONS

```
-i --input            DICOM filename / directory
-o --output           DICOM filename / directory
-r --recursive        recursively process (sub-)directories.
--continue            Do not stop when file found is not DICOM.
--root-uid            Root UID.
--resources-path      Resources path.
-k --key              Path to RSA Private Key.
-c --certificate       Path to Certificate.
```

### 4.4.3 encryption options

```
--des                DES.
--des3               Triple DES.
--aes128             AES 128.
--aes192             AES 192.
--aes256             AES 256.
```

### 4.4.4 dumb mode options

```
--empty %d,%d        DICOM tag(s) to empty
--remove %d,%d        DICOM tag(s) to remove
--replace %d,%d,%s    DICOM tag(s) to replace
```

### 4.4.5 general options

```
-h --help            print this help text and exit
-v --version         print version information and exit
-V --verbose         verbose mode (warning+error).
-W --warning         warning mode, print warning information
-E --error           error mode, print error information
-D --debug           debug mode, print debug information
```

#### 4.4.6 environment variable

```
GDCM_ROOT_UID Root UID
GDCM_RESOURCES_PATH path pointing to resources files (Part3.xml, ...)
```

### 4.5 Typical usage

#### 4.5.1 De-identification (anonymization, encrypt)

The only thing required for this operation is a certificate file (in PEM format).

```
$ gdcmanon --certificate certificate.pem -e original.dcm original_anonymized.dcm
```

You can use `-asn1` option from `gdcm dump` to dump the generated DataSet as ASN1 structure (see `gdcm dump(1)` for example).

#### 4.5.2 Re-identification (de-anonymization, decrypt)

The only thing required for this operation is a private key (in PEM format). It is required that the private key used for the re-identification process, was the actual private key used to generate the certificate file (certificate.pem) used during the de-identification step.

```
$ gdcmanon --key privatekey.pem -d original_anonymized.dcm original_copy.dcm
```

You can then check that `original.dcm` and `original_copy.dcm` are identical.

#### 4.5.3 Multiple files caveat

It is very important to understand the following section, when anonymizing more than one single file. When anonymizing multiple DICOM files, you are required to use the directory input. You cannot call multiple time the `gdcmanon` command line tool. Indeed the tool stores in memory during the process only a hash table of conversion so that each time a particular value is found it get always replaced by the same de-identified value (think: consistent Series Instance UID).

#### 4.5.4 Dumb mode

This functionality is not described in the DICOM standard. Users are advised that improper use of that mode is not recommended, meaning that important tag can be emptied/removed/replaced resulting in illegal/invalid DICOM file. Only use when you know what you are doing. If you delete a Type 1 attribute, chance is that your DICOM file will be not accepted in most DICOM third party viewer. Unfortunately this is often this mode that is implemented in popular DICOM Viewer, always prefer what the DICOM standard describes, and avoid the dumb mode.

The following example shows how to use dumb mode and achieve 5 operations at the same time:

- Empty the tag (0010,0010) Patient's Name,
- Empty the tag (0010,0020) Patient ID,
- Remove the tag (0010,0040) Patient's Sex
- Remove the tag (0010,1010) Patient's Age

- Replace the tag (0010,1030) Patient's Weight with the value '10'

You are required to check which DICOM attribute is Type 1 and Type 1C, before trying to '**Empty**' or '**Remove**' a particular DICOM attribute. For the same reason, you are required to check what are valid value in a replace operation.

```
$ gdcmanon --dumb --empty 10,10 --empty 10,20 --remove 10,40 --remove 10,1010 --replace 10,1030,10 012345.002.050
```

Multiple operation of `--dumb` mode can take place, just reuse the output of the previous operation. Always use `gdcmdump` on the input and output file to check what was actually achieved. You can use a diff program to check only what changed (see `gdcmdiff(1)` for example).

#### 4.5.4.1 Irreversible Anonymization

In some very rare cases, one would want to anonymize using the PS 3.15 mode so as to take benefit of the automatic conversion of all content that could contain Patient related information.

In the end all Patient related information has been removed and has been secretly stored in the 0400,0500 DICOM attribute. However to make sure that no-one ever try to break that security using brute-force algorithm, one want want to remove completely this DICOM attribute. This will make the DICOM:

- Completely free of any Patient related information (as per PS 3.15 specification)
- Remove any mean of people to brute force attack the file to find out the identity of the Patient

In this case one could simply do, as a first step execute the reversible anonymizer:

```
$ gdcmanon -c certificate.pem input.dcm anonymized_reversible.dcm
```

and now completely remove the DICOM attribute containing the secretly encrypted Patient related information:

```
$ gdcmanon --dumb --remove 400,500 --remove 12,62 --remove 12,63 anonymized_reversible.dcm anonymized_irreversible.dcm
```

#### Remarks

As mentionned in DICOM Sup 142, this anonymization is preferred over de-identification since: It is not required that the Encrypted Attributes Data Set be created; indeed, there may be circumstances where the Dataset is expected to be archived long enough that any contemporary encryption technology may be inadequate to provide long term protection against unauthorized recovery of identification

## 4.6 OpenSSL

On most system you can have access to OpenSSL to generate the Private Key/Certificate pair.

### 4.6.1 Generating a Private Key

Command line to generate a rsa key (512bit)

```
$ openssl genrsa -out CA_key.pem
```

Command line to generate a rsa key (2048bit)



```
$ openssl genrsa -out CA_key.pem 2048
```

Command line to generate a rsa key (2048bit) + passphrase

```
$ openssl genrsa -des3 -out CA_key.pem 2048
```

#### 4.6.2 Generating a Certificate

From your previously generated Private Key, you can now generate a certificate in PEM (DER format is currently not supported).

```
$ openssl req -new -key CA_key.pem -x509 -days 365 -out CA_cert.cer
```

## 4.7 DICOM Standard:

Page to the DICOM Standard:

<http://dicom.nema.org/>

The DICOM Standard at the time of releasing gdcmanon is:

<ftp://medical.nema.org/medical/dicom/2008/>

Direct link to PS 3.15-2008:

[ftp://medical.nema.org/medical/dicom/2008/08\\_15pu.pdf](ftp://medical.nema.org/medical/dicom/2008/08_15pu.pdf)

## 4.8 Warnings

Certain attributes may still contains Protected Health Information (PHI) after an anonymization step. This is typically the case for Patient's Address (0010,1040). The reason is that this particular attribute is not supposed to be in the composite IODs in the first place. DICOM Supp 142 includes it (however gdcmanon does not implement it).

## 4.9 SEE ALSO

**gdcconv(1), gdcmdump(1), gdcmdiff(1), openssl(1), dumpasn1(1)**

## 4.10 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 5

# Tool to convert DICOM to DICOM.

### 5.1 SYNOPSIS

```
gdcmconv [options] file-in file-out
```

### 5.2 DESCRIPTION

The **gdcmconv** command line program takes as input a DICOM file (file-in) and process it to generate an output DICOM file (file-out). The command line option dictate the type of operation(s) gdcmconv will use to generate the output file.

### 5.3 PARAMETERS

```
file-in    DICOM input filename
```

```
file-out    DICOM output filename
```

### 5.4 OPTIONS

#### 5.4.1 PARAMETERS

```
-i --input      DICOM filename
-o --output      DICOM filename
```

#### 5.4.2 OPTIONS

```
-X --explicit      Change Transfer Syntax to explicit.
-M --implicit      Change Transfer Syntax to implicit.
-U --use-dict       Use dict for VR (only public by default).
  --with-private-dict Use private dict for VR (advanced user only).
-C --check-meta     Check File Meta Information (advanced user only).
  --root-uid        Root UID.
  --remove-gl       Remove group length (deprecated in DICOM 2008).
  --remove-private-tags Remove private tags.
  --remove-retired  Remove retired tags.
```

### 5.4.3 image options

```
-l --apply-lut           Apply LUT (non-standard, advanced user only).
-P --photometric-interpretation %s Change Photometric Interpretation (when possible).
-w --raw                Decompress image.
-d --deflated            Compress using deflated (gzip).
-J --jpeg               Compress image in jpeg.
-K --j2k                Compress image in j2k.
-L --jpegls             Compress image in jpeg-ls.
-R --rle                Compress image in rle (lossless only).
-F --force              Force decompression/merging before recompression/splitting.
  --generate-icon        Generate icon.
  --icon-minmax %d,%d    Min/Max value for icon.
  --icon-auto-minmax     Automatically compute best Min/Max values for icon.
  --compress-icon        Decide whether icon follows main TransferSyntax or remains uncompressed.
  --planar-configuration [01] Change planar configuration.
-Y --lossy              Use the lossy (if possible) compressor.
-S --split %d           Write 2D image with multiple fragments (using max size)
```

### 5.4.4 JPEG options

```
-q --quality %*f        set quality.
```

### 5.4.5 JPEG-LS options

```
-e --lossy-error %*i    set error.
```

### 5.4.6 J2K options

```
-r --rate %*f           set rate.
-q --quality %*f        set quality.
-t --tile %d,%d         set tile size.
-n --number-resolution %d set number of resolution.
  --irreversible         set irreversible.
```

### 5.4.7 general options

```
-h --help               print this help text and exit
-v --version            print version information and exit
-V --verbose            verbose mode (warning+error).
-W --warning            warning mode, print warning information
-E --error              error mode, print error information
-D --debug              debug mode, print debug information
```

### 5.4.8 special options

```
-I --ignore-errors      convert even if file is corrupted (advanced users only, see disclaimers).
```

### 5.4.9 environment variable

```
GDCM_ROOT_UID Root UID
```

## 5.5 Simple usage

**gdcmmconv** is a great tool to convert broken DICOM implementation into properly parsable DICOM file. Usage is simply:

```
$ gdcmmconv input.dcm output.dcm
```

or if you prefer being explicit:

```
$ gdcmmconv -i input.dcm -o output.dcm
```

Even though **gdcmmconv** can overwrite directly on the same file (`input.dcm = output.dcm`), it is recommended that user should first convert into a different file to make sure the bug is properly handled by GDCM.

Typical cases where you would want to use **gdcmmconv** in its simple form:

- convert non-cp246 conforming file into conforming cp246,
- convert implicit little endian transfer syntax file meta header into proper explicit little endian transfer syntax,
- convert the GE-13 bytes bug,
- convert dual syntax file: implicit/explicit,
- convert Philips dual Little Endian/Big Endian file,
- convert GDCM 1.2.0 broken UN-2-bytes fields,
- &...
- All other broken files listed in the supported section.

When no option other is used, only the dataset is inspected. So encapsulated Pixel Data, for instance, is not inspected for well known bugs.

When doing this kind of work, this is usually a good idea to perform some kind of quality control, see **gdcmmconv** Quality Control section (down below).

## 5.6 Typical usage

### 5.6.1 File Meta Header

Running

```
$ gdcmmconv input.dcm output.dcm
```

Is not enough to recompute file meta header, when input file is buggy. You may want to use: `--check-meta`

```
$ gdcmmconv --check-meta input.dcm output.dcm
```

See typical cases such as: `GE_DLX-8-MONO2-PrivateSyntax.dcm` or `PICKER-16-MONO2-No_DicomV3_Preamble.dcm` from `gdcmmData`.

### 5.6.2 Conversion to Explicit Transfer Syntax

To convert a file that was written using Implicit Transfer Syntax into Explicit Transfer Syntax simply use:

```
$ gdcconv --explicit uncompressed.dcm compressed.dcm
```

### 5.6.3 Compressing to lossless JPEG

To compress an uncompressed DICOM file to a JPEG Lossless encapsulated format:

```
$ gdcconv --jpeg uncompressed.dcm compressed.dcm
```

### 5.6.4 Compressing to lossy JPEG

To compress an uncompressed DICOM file to a JPEG Lossy encapsulated format:

```
$ gdcconv --lossy --jpeg -q 90 uncompressed.dcm compressed.dcm
```

Note:

`-q` is just one of the many way to specify lossy quality, you need to inspect the other cmd line flag to specify

### 5.6.5 Compressing to lossless JPEG-LS

To compress an uncompressed DICOM file to a JPEG-LS Lossless encapsulated format:

```
$ gdcconv --jpegls uncompressed.dcm compressed.dcm
```

### 5.6.6 Compressing to lossy JPEG-LS

To compress an uncompressed DICOM file to a JPEG-LS Lossy encapsulated format:

```
$ gdcconv --lossy --jpegls -e 2 uncompressed.dcm lossy_compressed.dcm
```

Note:

`-e` (or `-lossy-error`) means that the maximum tolerate error is 2 for each pixel value

### 5.6.7 Compressing to lossless J2K

To compress an uncompressed DICOM file to a JPEG-2000 Lossless encapsulated format:

```
$ gdcconv --j2k uncompressed.dcm compressed.dcm
```

### 5.6.8 Compressing to lossy J2K

To compress an uncompressed DICOM file to a JPEG-2000 Lossy encapsulated format:

```
$ gdcconv --lossy -q 55,50,45 --j2k uncompressed.dcm lossy_compressed.dcm
```

Note:

`-q` is just one of the many way to specify lossy quality, you need to inspect the other cmd line flag to specify

### 5.6.9 Compressing to lossless RLE

To compress an uncompressed DICOM file to a RLE Lossless encapsulated format:

```
$ gdcconv --rle uncompressed.dcm compressed.dcm
```

There is no such thing as lossy RLE compression.

### 5.6.10 Split encapsulated DICOM:

To split an encapsulated stream into smaller chunk (1024 bytes each):

```
$ gdcconv --split 1024 rle.dcm rle_1024.dcm
```

If an odd number of bytes is passed it will be rounded down to the next even number (eg. 1025 -> 1024) since DICOM only allow even number for Value Length.

### 5.6.11 Forcing (re)compression

Sometime it is necessary to use the `-force` option. By default when user specify `-j2k` and input file is already in JPEG 2000 encapsulated DICOM format then no operation takes places. By using `-force` you make sure that (re)compression operation takes places.

Real life example of why you would use `-force`:

- When Pixel Data is missing data / is padded with junk
- When you would like to make sure GDCM can handle decompression & recompression cycle

### 5.6.12 Decompressing a Compressed DICOM

```
$ gdcconv --raw compressed.dcm uncompressed.dcm
```

### 5.6.13 Compressing an uncompressed Icon

By default when compressing a DICOM Image file, `gdcconv` will not compress the icon. A user option needs to be turned on to explicitly force the compression of the Icon Image Sequence Pixel Data

For example, by default we will not compress the Icon Image Sequence Pixel Data attribute:

```
$ gdcconv --jpeg gdcData/simpleImageWithIcon.dcm uncompressed_icon.dcm
```

In the following example we will explicitly compress the Icon Image Sequence Pixel Data attribute. In that case the same Transfer Syntax is being used for both the main Pixel Data and the Pixel Data from the Icon Image Sequence:

```
$ gdcconv --jpeg --compress-icon gdcData/simpleImageWithIcon.dcm compressed_icon.dcm
```

### 5.6.14 Generating an Icon

For some application it might be necessary to produce a small preview of the main image to be able to quickly load that short preview instead of the main image. In that case:

```
$ gdcconv --raw --generate-icon gdcData/test.acr test_icon.dcm
```

In some cases the main Pixel Data element is expressed as pixel defined on 16bits. Since Icon can only store at most pixel of size 8bits, a rescale operation needs to take place. In order to properly select a better interval for doing the rescale operation user can specify the min max used for the rescale operation:

```
$ gdcconv --raw --generate-icon --icon-minmax 0,192 gdcData/012345.002.050.dcm icon_minmax.dcm
```

### 5.6.15 Changing the planar Configuration

Often RLE files are compressed using a different Planar Configuration (RRR ... GGG... BBB...) instead of the usual triplet (RGB ... RGB ... RGB ). So upon decompression the Planar Configuration is 1. This is perfectly legal in DICOM, however this is unconventional, and thus it may be a good idea to also change the planar configuration and set it to the default :

```
$ gdcconv --raw --planar-configuration 0 compressed.dcm uncompressed1.dcm
```

To reinvert the planar configuration of file 'uncompressed1.dcm', simply do:

```
$ gdcconv --raw --planar-configuration 1 uncompressed1.dcm uncompressed2.dcm
```

## 5.7 Lossless Conversion

When talking about lossless conversion, there is an ambiguity that need to be understood. To achieve higher compression ratio, the RGB color space is usually not used, in favor of a YBR one. Changing from one color space to the other is (bit level) not lossless.

For more detail, see what are the true lossless transformations as described:

[http://gdc.sourceforge.net/wiki/index.php/Color\\_Space\\_Transformations](http://gdc.sourceforge.net/wiki/index.php/Color_Space_Transformations)

## 5.8 Quality Control

One important part when using gdcconv it to have a way to quality control the output.

You can use 3rd party tool to check the output of gdcconv is correct.

### 5.8.1 DCMTK / dicom3tools

Using another DICOM implementation such as the one from DCMTK or dicom3tools can be a good process to check the output of gdcconv.

- For DCMTK use: dcmdump
- For dicom3tools use: dcdump

For reference, gdcconv --raw will act as dcmdjpeg +cn +px, since it never tries to convert color space.



### 5.8.2 VIM: vimdiff

You can setup your favorite editor to compare the output, for instance in vim:

```
autocmd BufReadPre *.dcm set ro
autocmd BufReadPost *.dcm silent %!gdcm dump -M +uc "%"
```

then simply do:

```
$ vimdiff input.dcm output.dcm
```

### 5.8.3 vbindiff

On UNIX you can visually compare binary file using the vbindiff command:

```
$ vbindiff input.dcm output.dcm
```

## 5.9 SEE ALSO

**gdcmdump(1)**, **gdcmmraw(1)**, **gdcminfo(1)**, **gdcmdiff(1)**

## 5.10 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 6

# dumps differences of two DICOM files

### 6.1 SYNOPSIS

```
gdcmdiff [options] file1 file2
```

### 6.2 DESCRIPTION

The **gdcmdiff** command line program takes as input two DICOM files: file1 and file2.

### 6.3 PARAMETERS

```
file1    DICOM input filename  
file2    DICOM output filename
```

### 6.4 OPTIONS

#### 6.4.1 OPTIONS

```
-m      --meta          Compare metainformation. Default is off.  
-t <n>  --truncate <n> String values trimmed to n characters.
```

#### 6.4.2 general options

```
-h      --help          print this help text and exit  
  
-v      --version       print version information and exit  
  
-V      --verbose       verbose mode (warning+error).  
  
-W      --warning       warning mode, print warning information
```

```
-E  --error  
    error mode, print error information  
  
-D  --debug  
    debug mode, print debug information
```

## 6.5 Simple usage

**gdcmdiff** is a great tool to produce a diff in between two DICOM files. Usage is simply:

```
$ gdcmdiff input1.dcm input2.dcm
```

## 6.6 SEE ALSO

**gdcmdump(1)**, **gdcminfo(1)**

## 6.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 7

**dumps a DICOM file, it will display the structure and values contained in the specified DICOM file.**

### 7.1 SYNOPSIS

```
gdcmdump [options] dcm_file
gdcmdump [options] dcm_directory
```

### 7.2 DESCRIPTION

The **gdcmdump** command line program dumps a DICOM file to the console. For those familiar with dcmdump (DCMTK) output, gdcmdump has some minor differences. Namely:

- For Implicit Transfer Syntax gdcmdump will print ?? instead of the dictionary VR

gdcmdump has a limited private dictionary that is used to lookup private element whenever possible.

### 7.3 PARAMETERS

```
dcm_file          DICOM input filename
dcm_directory     DICOM input directory
```

### 7.4 OPTIONS

#### 7.4.1 OPTIONS

-x --xml-dict	generate the XML dict (only private elements for now).
-r --recursive	recursive (input is a directory)
-d --dump	dump value (limited use).
-p --print	print value instead of simply dumping (default).
-c --color	print in color.
-C --csa	print SIEMENS CSA Header (0029,[12]0,SIEMENS CSA HEADER).
-P --pdb	print GEMS Protocol Data Block (0025,1b,GEMS_SERS_01).
--elscint	print ELSCINT Protocol Information (01f7,26,ELSCINT1).
--vepro	print VEPRO Protocol Information (0055,20,VEPRO VIF 3.0 DATA).

```

                or VEPRO Protocol Information (0055,20,VEPRO VIM 5.0 DATA).
--sds           print Philips MR Series Data Storage (1.3.46.670589.11.0.0.12.2) Information (2005,32,Philips)
-A --asn1       print encapsulated ASN1 structure >(0400,0520).
--map-uid-names map UID to names.

```

## 7.4.2 general options

```

-h --help
    print this help text and exit

-v --version
    print version information and exit

-V --verbose
    verbose mode (warning+error).

-W --warning
    warning mode, print warning information

-E --error
    error mode, print error information

-D --debug
    debug mode, print debug information

```

## 7.4.3 special options

```

-I --ignore-errors  dumps even if file is corrupted (advanced users only, see disclaimers).

```

## 7.5 Typical usage

### 7.5.1 Printing Implicit Transfer Syntax

The VR are not found in the file, thus are presented with a "(??)", and right next to it (if found) the correct VR.

Eg.:

```
$ gdcmdump GE_DLX-8-MONO2-PrivateSyntax.dcm
```

```

# Dicom-File-Format
\&...
(0008,0000) ?? (UL) 434                                # 4,1 Generic Group Length
(0008,0005) ?? (CS) [ISO_IR 100]                        # 10,1-n Specific Character Set
(0008,0008) ?? (CS) [ORIGINAL\\PRIMARY\\SINGLE PLANE ]   # 30,2-n Image Type
(0008,0016) ?? (UI) [1.2.840.10008.5.1.4.1.1.12.1]      # 28,1 SOP Class UID
(0008,0018) ?? (UI) [1.2.840.113619.2.16.1.0.906539207.1.24207] # 42,1 SOP Instance UID
(0008,0020) ?? (DA) [19980923]                          # 8,1 Study Date
(0008,0021) ?? (DA) [19980923]                          # 8,1 Series Date
(0008,0022) ?? (DA) [19980923]                          # 8,1 Acquisition Date
(0008,0023) ?? (DA) [19980923]                          # 8,1 Content Date
(0008,0030) ?? (TM) [101229.000]                        # 10,1 Study Time
(0008,0031) ?? (TM) [101229.000]                        # 10,1 Series Time
(0008,0032) ?? (TM) [102653.000]                        # 10,1 Acquisition Time
(0008,0033) ?? (TM) [102653.000]                        # 10,1 Content Time
\&...

```

### 7.5.2 Print Private Attributes

GDCM has a limited private dictionary. Whenever possible, it will try to lookup the private data element.

```
$ gdcmdump 012345.002.050.dcm
```

```
\&...
(0009,0010) LO [GEMS_IDEN_01] # 12,1 Private Creator
(0009,1001) LO [GE_GENESIS_FF ] # 14,1 Full fidelity
(0009,1002) SH [MRCV] # 4,1 Suite id
(0009,1004) SH [SIGNA ] # 6,1 Product id
(0009,1027) SL 985968524 # 4,1 Image actual date
(0009,1030) SH [19356UMR2 ] # 10,1 Service id
(0009,1031) SH [999 ] # 4,1 Mobile location number
(0009,10e3) UI [1.2.840.113619.1.1.4.1762386977] # 32,1 Equipment UID
(0009,10e6) SH [08] # 2,1 Genesis Version - now
(0009,10e7) UL 2757786872 # 4,1 Exam Record checksum
(0009,10e9) SL 985968523 # 4,1 Actual series data time stamp
\&...
(0019,0000) UL 1208 # 4,1 Generic Group Length
(0019,0010) LO [GEMS_ACQU_01] # 12,1 Private Creator
(0019,100f) DS [424.399994] # 10,1 Horiz. Frame of ref.
(0019,1011) SS 0 # 2,1 Series contrast
\&...
(0019,10e0) DS [0.000000] # 8,1 User data 24 {# DTI Diffusion Dir., relea
(0019,10e2) DS [0.000000] # 8,1 Velocity Encode Scale
(0019,10f2) SS 0 # 2,1 Fast phases
(0019,10f9) DS [98] # 2,1 Transmit gain
\&...
(0021,0000) UL 372 # 4,1 Generic Group Length
(0021,0010) LO [GEMS_RELA_01] # 12,1 Private Creator
(0021,1003) SS 0 # 2,1 Series from which Prescribed
\&...
```

### 7.5.3 SIEMENS CSA Header

Using this option it is possible to dump as a readable text what is contained in the private attribute as found in typical SIEMENS MR DICOM file.

Eg.:

```
$ gdcmdump --csa MR_SIEMENS_forceLoad29-1010_29-1020.dcm
```

```
(0029,0010)siemens csa header
Image shadow data (0029,xx10)

0 - 'EchoLinePosition' VM 1, VR IS, SyngoDT 6, NoOfItems 6, Data '64      '
1 - 'EchoColumnPosition' VM 1, VR IS, SyngoDT 6, NoOfItems 6, Data '64      '
2 - 'EchoPartitionPosition' VM 1, VR IS, SyngoDT 6, NoOfItems 6, Data '32      '
3 - 'UsedChannelMask' VM 1, VR UL, SyngoDT 9, NoOfItems 6, Data '255      '
4 - 'Actual3DImaPartNumber' VM 1, VR IS, SyngoDT 6, NoOfItems 0, Data
5 - 'ICE_Dims' VM 1, VR LO, SyngoDT 19, NoOfItems 6, Data 'X_1_1_1_1_1_31_1_1_1_1_19'
6 - 'B_value' VM 1, VR IS, SyngoDT 6, NoOfItems 6, Data '0      '
7 - 'Filter1' VM 1, VR IS, SyngoDT 6, NoOfItems 0, Data
8 - 'Filter2' VM 1, VR IS, SyngoDT 6, NoOfItems 0, Data
\&...
```

### 7.5.4 GEMS Protocol Data Block

Using this option it is possible to dump as a readable text what is contained in the private attribute as found in typical GEMS MR DICOM file.

Protocol Data Block : 0025,xx1b,GEMS\_SERS\_01

```
$ gdcmdump --pdb GE_MR_0025xx1bProtocolDataBlock.dcm
```

```
ENTRY "Head First"
POSITION "Supine"
ANREF "NA"
COIL "HEAD"
PLANE "OBLIQUE"
SEDESCFLAG "1"
SEDESC "AX FSE T2"
IMODE "2D"
PSEQ "FSE-XL"
IOPT "FC, EDR, TRF, Fast"
PLUG "22"
FILTCHOICE "None"
BWRT "-1"
TRICKSIMG "1"
TAG_SPACE "7"
TAG_TYPE "None"
\&...
```

### 7.5.5 ELSCINT Protocol Information

Using this option it is possible to dump as a readable text what is contained in the private attribute as found in typical ELSCINT CT DICOM file.

ELSCINT Protocol Information: (01f7,26,ELSCINT1)

```
$ gdcmdump --elscint ELSCINT1_ProtocolInformation.dcm
```

```
ELSCINT1 Dumping info from tag (01f7,26,elscint1)
```

```
ELSCINT1/Item name: []
  ApprovedStep [yes]
  RefSurview [1\0]
  STD-first-img-pos [11.5]
  current-step [yes]
  ntimed-steps [0]
  orig-n-slices [390]
  protocol-file [Head_Multi_1032_usr.proc]
  protocol-name [FACE-TRAUMA/Head/Hx]
  protocol-path [/usr/diamond.root/spr/]
  protocol-step [1]
  protocol-version [2.51]
```

```
ELSCINT1/Item name: [doseright]
```

```
  ACS [n/a]
  ACS-bed-position [0]
  ACS-calc-mas [0]
  ACS-ig-parameter [0]
  ACS-learn-allowed [no]
  ACS-water-radius [-1.000000]
  ACS-water-radius-scan [-1]
\&...
```

### 7.5.6 VEPRO Protocol Information

Using this option it is possible to dump as a readable text what is contained in the private attribute as found in typical VEPRO CT DICOM file.

ELSCINT Protocol Information: (0055,20,VEPRO VIM 5.0 DATA)



```
$ gdcmdump --vepro VEPRO_ProtocolInformation.dcm

VIMDATA2: (0055,20,VEPRO VIM 5.0 DATA)
  ID: VIM
  Version: 5.0
  UserName:
  UserAdress1: Name of Institution
  UserAdress2: Street of Institution
  UserAdress3: City of Institution
  UserAdress4:
  UserAdress5:
  RecDate: 20101001
  RecTime: 211321
  RecPlace:
  RecSource: DICOM Distributor
  DF1: P-09/10-41808
  DF2: Sultana Razia
  DF3: 19411001
  DF4: F
  DF5:
  DF6:
  DF7:
  DF8: CT Scan Brain without Contrast
  DF9: 10/10-0034873
  DF10: 10/10-00348
  DF11:
  DF12:
  DF13:
  DF14: Head 0.5
  DF15: 4
  DF16:
  DF17:
  DF18:
  DF19:
  DF20:
  StudyUID: 1.2.392.200036.9116.2.6.1.48.1214228007.1285934880.206831
  SeriesUID: 1.2.392.200036.9116.2.6.1.48.1214228007.1285935201.938653
  Modality: CT
```

### 7.5.7 Philips Private MR Series Data Storage (1.3.46.670589.11.0.0.12.2)

Using this option it is possible to dump as a readable text what is contained in the private attribute as found in typical Philips Private MR Series Data Storage file.

PMS Series Data Storage (2005,32,Philips MR Imaging DD 002)

```
$ gdcmdump --sds PMS_SeriesDataStorage.dcm

\&...
PMS/Item name: [PDF_CONTROL_GEN_PARS/IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_CONTROL_PREP_PARS /IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_CONTROL_RECON_PARS/IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_CONTROL_SCAN_PARS /IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_EXAM_PARS /IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_HARDWARE_PARS /IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_PREP_PARS /IEEE_PDF/Y ]
\&...
PMS/Item name: [PDF_SPT_PARS/IEEE_PDF/Y ]
```

```

SP_scan_resol [256\256] # 2
SP_pda_profiles [0\0] # 2
SP_filter [324074] # 1
SP_analyse_with_iqt [0] # 1
SP_main_system_type [3] # 1
SP_gradient_system [6] # 1
SP_coil_type [2\2\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_coil_id [2\34\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_coil_part [0\0\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_act_q [0\0\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_act_coil_freq [0\0\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_coil_m_pos [255\255\255\0\0\0\0\0\0\0\0\0\0\0\0\255] # 16
SP_coil_t_pos [255\128\255\0\0\0\0\0\0\0\0\0\0\0\0\255] # 16
SP_surface_coil_con [0\1\0\0\0\0\0\0\0\0\0\0\0\0\0\0] # 16
SP_proton_freq [127801349] # 1
SP_tm_result [2\2\2\2\2\2\2\2\2\2\2\2\2\2\2\2] # 16
SP_f0_result [0] # 1
SP_as_result [0] # 1
SP_po_result [0] # 1
SP_rg_result [0] # 1
SP_dc_result [0] # 1
SP_ph_result [0] # 1
\&...

```

## 7.5.8 Encapsulated ASN1 Structure

This option is mainly used for dumping the ASN1 structure of the encrypted Attribute (0040,0520)

```
$ gdcmdump encrypted.dcm
```

```

\&...
(0400,0500) SQ # u/1,1 Encrypted Attributes Sequence
  (fffe,e000) na (Item with undefined length)
    (0400,0510) UI [1.2.840.10008.1.2] # 18,1 Encrypted Content Transfer Syntax UID
    (0400,0520) OB 30\82\03\ba\06\09\2a\86\48\55\04\08\13 # 958,1 Encrypted Content
  (fffe,e00d)
(fffe,e0dd)
\&...

```

```
$ gdcmdump --asn1 encrypted.dcm
```

```

0:d=0 hl=4 l= 954 cons: SEQUENCE
4:d=1 hl=2 l= 9 prim: OBJECT :pkcs7-envelopedData
15:d=1 hl=4 l= 939 cons: cont [ 0 ]
19:d=2 hl=4 l= 935 cons: SEQUENCE
23:d=3 hl=2 l= 1 prim: INTEGER :00
26:d=3 hl=4 l= 366 cons: SET
30:d=4 hl=4 l= 362 cons: SEQUENCE
34:d=5 hl=2 l= 1 prim: INTEGER :00
37:d=5 hl=2 l= 82 cons: SEQUENCE
39:d=6 hl=2 l= 69 cons: SEQUENCE
41:d=7 hl=2 l= 11 cons: SET
43:d=8 hl=2 l= 9 cons: SEQUENCE
45:d=9 hl=2 l= 3 prim: OBJECT :countryName
50:d=9 hl=2 l= 2 prim: PRINTABLESTRING :AU
54:d=7 hl=2 l= 19 cons: SET
56:d=8 hl=2 l= 17 cons: SEQUENCE
58:d=9 hl=2 l= 3 prim: OBJECT :stateOrProvinceName
63:d=9 hl=2 l= 10 prim: PRINTABLESTRING :Some-State
75:d=7 hl=2 l= 33 cons: SET
77:d=8 hl=2 l= 31 cons: SEQUENCE
79:d=9 hl=2 l= 3 prim: OBJECT :organizationName
84:d=9 hl=2 l= 24 prim: PRINTABLESTRING :Internet Widgits Pty Ltd
110:d=6 hl=2 l= 9 prim: INTEGER :AC966D88787A51B4

```

```

121:d=5 hl=2 l= 13 cons: SEQUENCE
123:d=6 hl=2 l= 9 prim: OBJECT :rsaEncryption
134:d=6 hl=2 l= 0 prim: NULL
136:d=5 hl=4 l= 256 prim: OCTET STRING [HEX DUMP]:822368070285AD756C962ECB973514B291F946...
396:d=3 hl=4 l= 558 cons: SEQUENCE
400:d=4 hl=2 l= 9 prim: OBJECT :pkcs7-data
411:d=4 hl=2 l= 29 cons: SEQUENCE
413:d=5 hl=2 l= 9 prim: OBJECT :aes-256-cbc
424:d=5 hl=2 l= 16 prim: OCTET STRING [HEX DUMP]:3B49AFE71749F2BFF1519EBAEA95A393
442:d=4 hl=4 l= 512 prim: cont [ 0 ]

```

## 7.6 SEE ALSO

**gdcmdump(1)**, **gdcmrw(1)**, **gdcmanon(1)**

## 7.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 8

# Tool to generate a DICOMDIR file from a File-Set.

### 8.1 SYNOPSIS

```
gdcmgendir [options] file-in file-out
```

### 8.2 DESCRIPTION

### 8.3 PARAMETERS

file-in    DICOM input filename

file-out   DICOM output filename

### 8.4 OPTIONS

#### 8.4.1 Parameters

#### 8.4.2 OPTIONS

-i --input	DICOM filename or directory
-o --output	DICOM filename or directory
-r --recursive	recursive.
--descriptor	descriptor.
--root-uid	Root UID.

#### 8.4.3 general options

-h    --help	print this help text and exit
-v    --version	print version information and exit

```
-V  --verbose
    verbose mode (warning+error).

-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

#### 8.4.4 environment variable

```
GDCM_ROOT_UID Root UID
```

### 8.5 Typical usage

### 8.6 NOTE

One may have to run some preliminary steps in order to get gdcmgendir to generate the DICOMDIR file. Namely two steps:

- Batch renaming of the DICOM filename into something compatible with ISO 9660 filename convention
- Convert all DICOM file into the Explicit VR Little Endian Uncompressed (1.2.840.10008.1.2.1)

Step 1. can be solved in a numerous way. Eg. on UNIX environment this could either be solved using the `mkisofs` command line tool. Filenames should not contains any extension since the VR CS does not allow for the '.' character. Only upper case, digit 0-9, the space ' ' and the underscore '\_' character are valid in VR CS, with a maximum of 8 bytes. Another simple tool that can be handy is 'rename' in conjunction with 'basename'.

Step 2. can simply be achieved using the `gdcconv` command line tool:

```
$ for i in `ls IMG*`; do gdcconv --raw --force $i /tmp/out/$i; done
```

### 8.7 SEE ALSO

**gdcconv(1)**, **gdcmanon(1)**, **rename(1)**, **mkisofs(1)**

### 8.8 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 9

# Manipulate DICOM image file.

`gdcmimg` is a low level tool to allow de-/encapsulation from/to DICOM image. This tool does not understand Transfer Syntax conversion. It will encapsulate the raw data as-is. This has some impact in some cases, see special warnings below.

It is important to note that `gdcmimg` can only encapsulate proper input file, for instance JPG and or JP2 are accepted since an associated DICOM Transfer Syntax can be found. However input such as TIFF and/or PNG are not, since DICOM does not support those. See instead a tool such as `gdcm2vtk`.

### 9.1 SYNOPSIS

```
gdcmimg [options] file-in file-out
```

### 9.2 DESCRIPTION

The **`gdcmimg`** command line tool can be used in two fashions:

- 1. Converting a recognized file format into its encapsulated DICOM counterpart,
- 2. Anonymizing a rectangular portion of a DICOM file.

### 9.3 PARAMETERS

```
file-in    input filename (non-DICOM)
```

```
file-out   DICOM output filename
```

### 9.4 OPTIONS

#### 9.4.1 PARAMETERS

```
-i --input      Input filename  
-o --output     Output filename
```

## 9.4.2 OPTIONS

```

--endian %s      Endianness (LSB/MSB).
-d --depth %d    Depth (Either 8/16/32 or BitsAllocated eg. 12 when known).
--sign %s        Pixel sign (0/1).
--spp %d         Sample Per Pixel (1/3).
-s --size %d,%d  Size.
-C --sop-class-uid SOP Class UID (name or value).
-T --study-uid   Study UID.
-S --series-uid  Series UID.
--root-uid       Root UID.

```

## 9.4.3 fill options

```

-R --region %d,%d Region.
-F --fill %d      Fill with pixel value specified.

```

## 9.4.4 general options

```

-h --help
    print this help text and exit

-v --version
    print version information and exit

-V --verbose
    verbose mode (warning+error).

-W --warning
    warning mode, print warning information

-E --error
    error mode, print error information

-D --debug
    debug mode, print debug information

```

## 9.4.5 environment variable

```
GDCM_ROOT_UID Root UID
```

## 9.5 Supported File Format (appropriate file extension) gdcming

will base it's conversion process based on the file extension. Follows the list of recognized file extension. When no extension is found, DICOM file is assumed.

input format

```

* RAW      (raw, rawl, gray, rgb)
* RLE      (rle)
* PNM      (pgm, pnm, ppm)
* JPEG-LS  (jls)
* JPEG 2000 (jp2, j2k, j2c, jpx, jpc)
* JPEG     (jpg, jpeg, ljpg, ljpeg)
* DICOM    ()

```

output format:



```
* PGM      (pgm, pnm, ppm)
* DICOM    ()
```

For RAW file format, you should take special care of the `--endian` option. For the (old) JPEG file format, both the lossy and lossless format are supported, user should pay attention to the `--sign` option. For file format such as RLE or RAW, user is expected to fill in information required to find the dimension and type of input data as there is no other way to find this information. For all other file format, the properties are derived from the file format itself.

PNM file are supposed to be big endian (important for depth > 8)

## 9.6 Typical usage

### 9.6.1 Remove a rectangular part of the image

To fill the region [0,100]x[0,100] of a DICOM image simply do:

```
$ gdcimg --fill 0 --region 0,100,0,100 -i input.dcm -o output_black.dcm
```

Warning: if the Pixel Data is compressed, the image is first decompressed so that pixel can be set to 0, but it is not re-compressed.

### 9.6.2 Convert RAW to DICOM

Recognized extension is `.raw`, `.rawl`, `.gray` or `.rgb` (case insensitive)

```
$ gdcimg --size 512,512 --depth 16 -i input.raw -o output.dcm
```

the image will be a Secondary Capture.

When the input is 3 component, one need to specify explicitly the Samples Per Pixel:

```
$ gdcimg --size 512,512 --spp 3 input_rgb.raw output_rgb.dcm
```

When the filename contains `.rgb` as file extension output is automatically recognized as RGB no need to specify `--spp`

```
$ gdcimg --size 512,512 input.rgb output_rgb.dcm
```

You can use the `dd` cmd line to skip any header you would like to discard, for instance, if you would like to skip the first 108 bytes, simply do:

```
$ dd skip=108 bs=1 if=input.raw of=output.raw
```

`.raw` and `.rawl` extension are equivalent. You need to explicitly specify the endianness manually:

```
$ gdcimg --endian MSB --size 512,512 --depth 16 -i input.raw -o output.dcm
```

or

```
$ gdcimg --endian LSB --size 512,512 --depth 16 -i input.raw -o output.dcm
```

### 9.6.3 Convert PGM/PNM/PPM to DICOM

Recognized extensions are .pgm, .pnm, .ppm (case insensitive)

```
$ gdcming -i input.pgm -o output.dcm
```

the image will be a Secondary Capture

### 9.6.4 Convert RLE to DICOM

Recognized extension is .rle (case insensitive)

```
$ gdcming --size 512,512 --depth 16 -i input.rle -o output.dcm
```

the image will be a Secondary Capture

### 9.6.5 Convert JPEG to DICOM

Recognized extensions are .jpg, .jpeg, .ljpg, .ljpeg (case insensitive)

```
$ gdcming -i input.ljpeg -o output.dcm
```

the image will be a Secondary Capture

### 9.6.6 Convert J2K to DICOM

Recognized extensions are .j2k, .jp2, .jpc, .jpx, .j2c (case insensitive)

```
$ gdcming -i input.j2k -o output.dcm
```

the image will be a Secondary Capture.

All Pixel informations (Bits Stored/Allocated...) will be derived from the image itself, and not from the command line options.

### 9.6.7 Specifying a SOP Class UID

Instead of the default Secondary Capture Image Storage, one may want to specify, say VL Photographic Image Storage.

```
$ gdcming --sop-class-uid 1.2.840.10008.5.1.4.1.1.77.1.4 input.jpg output.dcm
```

## 9.7 Multiple Files

gdcming handle nicely a set of files (for instance jpeg):

```
$ gdcming 1.jpg 2.jpg 3.jpg 4.jpg output.dcm
```

## 9.8 Start Offset

In some case, one may want to create a 2D slice from an arbitrary volume (e.g 3D). In which case `--offset` becomes handy:

```
$ gdcming --offset 4954104330 --size 1673,1673 Input3D_1673_1673_1775.raw slice_1770.dcm
```

## 9.9 Warning

There are a couple of issues with `gdcming` implementation:

For RAW file, one should pay attention that when using `--endian MSB` the Pixel Data will be encapsulated as is (not touched by `gdcming`). Therefore the only possible transfer syntax available is Implicit VR Big Endian DLX (G.E Private). GDCM does handle this private Transfer Syntax. So if you need to convert this Transfer Syntax to another one (and allow Pixel Data manipulation), you can use:

```
$ gdcconv --raw --force input_big_endian_dlx.raw -o output_implicit_vr_little_endian.dcm
```

For JFIF file and JP2 file (with header) the header is copied into the Pixel Data element which is illegal for JP2. Use `gdcconv` to properly re-encode a JP2/JFIF file into J2K/JPG.

```
$ gdcming input.jp2 output_jp2.dcm
$ gdcconv --j2k --force output_jp2.dcm output_j2k.dcm
```

For RLE file, no check is done for crossing the row boundary. It is recommended to use `gdcconv -rle` to re-encode into a proper RLE file in case of doubt.

Of course if the compression is not ok with your setup, you can always de-encapsulated the DICOM file (typically JPEG) to a non-encapsulated form, using `gdcconv`:

```
$ gdcconv --raw input_jpeg.dcm output_raw.dcm
```

## 9.10 SEE ALSO

**`gdcmdump(1)`, `gdc2vtk(1)`, `gdc2raw(1)`, `convert(1)`, `dd(1)`**

## 9.11 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 10

# Display meta info about the input DICOM file.

### 10.1 SYNOPSIS

```
gdcminfo [options] file-in
```

### 10.2 DESCRIPTION

The **gdcminfo** command line program takes as input a DICOM file, or a directory and process it to extract meta-information about the DICOM file processed.

### 10.3 PARAMETERS

```
file-in    DICOM input filename
```

### 10.4 OPTIONS

#### 10.4.1 OPTIONS

-r --recursive	recursive.
-d --check-deflated	check if file is proper deflated syntax.
--resources-path	Resources path.
--md5sum	Compute md5sum of Pixel Data attribute value.
--check-compression	check the encapsulated stream compression (lossless/lossy).

#### 10.4.2 general options

-h	--help	print this help text and exit
-v	--version	print version information and exit
-V	--verbose	verbose mode (warning+error).

```
-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

### 10.4.3 environment variable

GDCM\_RESOURCES\_PATH path pointing to resources files (Part3.xml, ...)

## 10.5 Simple usage

### 10.5.1 gdcmdata

Using data from gdcmdata:

```
$ gdcminfo gdcmdata/012345.002.050.dcm
```

```
MediaStorage is 1.2.840.10008.5.1.4.1.1.4 [MR Image Storage]
NumberOfDimensions: 2
Dimensions: (256,256)
Origin: (-85,21.6,108.7)
Spacing: (0.664062,0.664062,1.5)
DirectionCosines: (1,0,0,0,0,-1)
Rescale Intercept/Slope: (0,1)
SamplesPerPixel      :1
BitsAllocated        :16
BitsStored           :16
HighBit              :15
PixelRepresentation:1
Orientation Label: CORONAL
```

### 10.5.2 Davie Clunie datasets:

Using data from David Clunie datasets:

```
$ gdcminfo BRTUM001.dcm
```

```
MediaStorage is 1.2.840.10008.5.1.4.1.1.4.1 [Enhanced MR Image Storage]
NumberOfDimensions: 3
Dimensions: (256,256,15)
Origin: (40,-105,105)
Spacing: (0.820312,0.820312,6)
DirectionCosines: (0,1,0,0,0,-1)
Rescale Intercept/Slope: (0,1)
SamplesPerPixel      :1
BitsAllocated        :16
BitsStored           :16
HighBit              :15
PixelRepresentation:1
Orientation Label: SAGITTAL
```

### 10.5.3 Checking the md5sum of the Pixel Data

After compressing a DICOM file (see `gdcmconv`) using a lossless compression algorithm, it is fairly easy to compare the two files for differences at DICOM attribute level. However one operation is slightly easier to do: how to make sure the compression was actually lossless ? In this case one could use the `--md5sum` operation.

Take an uncompressed DICOM image file:

```
$ gdcminfo --md5sum SIEMENS_ImageLocationUN.dcm
```

The tool return: 0621954acd5815e0b4f7b65fcc6506b1

Now compress this file:

```
$ gdcmmconv --jpegls SIEMENS_ImageLocationUN.dcm lossless_compressed.dcm
```

and then check again the md5sum:

```
$ gdcminfo --md5sum lossless_compressed.dcm
```

The tool return: 0621954acd5815e0b4f7b65fcc6506b1

### 10.5.4 Checking if Pixel Data is lossless

In some environment one wish to check whether or not the DICOM file is lossless. It is fairly easy to do that in most cases. Only in two occasion this is not clear from the sole DICOM Attribute. When the Transfer Syntax is JPEG 2000 Image Compression (1.2.840.10008.1.2.4.91) and when the Transfer Syntax is JPEG-LS Lossy (Near-Lossless) Image Compression (1.2.840.10008.1.2.4.81).

In this case, the only solution is to open the Pixel Data element, read the specific JPEG header and check whether or not the JPEG transformation was lossless or not:

```
$ gdcminfo --check-compression gdcmmData/MAROTECH_CT_JP2Lossy.dcm
```

The tool returns: "Encapsulated Stream was found to be: lossy"

## 10.6 SEE ALSO

`gdcmdump(1)`, `gdcmraw(1)`, `gdcmconv(1)`

## 10.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre





## Chapter 11

# Tool to convert PAPYRUS 3.0 to DICOM.

### 11.1 SYNOPSIS

```
gdcmconv [options] file-in file-out
```

### 11.2 DESCRIPTION

The **gdcmconv** command line program takes as input a PAPYRUS 3.0 file (file-in) and process it to generate an output (pseudo) DICOM file (file-out). The command line option dictate the type of operation(s) gdcmconv will use to generate the output file.

### 11.3 PARAMETERS

```
file-in    DICOM input filename
```

```
file-out    DICOM output filename
```

### 11.4 OPTIONS

#### 11.4.1 PARAMETERS

```
-i --input      DICOM filename
-o --output     DICOM filename
```

#### 11.4.2 OPTIONS

```
-S --split      Split multiframe PAPYRUS 3.0 into multiples DICOM files
--decomp-pap3   Use PAPYRUS 3.0 for decompressing (can be combined with --split).
--check-iop     Check that the Image Orientation (Patient) Attribute is ok (see --split).
```

#### 11.4.3 general options

```
-h --help
```

```
    print this help text and exit

-v  --version
    print version information and exit

-V  --verbose
    verbose mode (warning+error).

-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

#### 11.4.4 environment variable

GDCM\_ROOT\_UID Root UID

### 11.5 Simple usage

**gdcmap3** is a great tool to convert broken PAPYRUS 3.0 implementation into properly parsable DICOM file. Usage is simply:

```
$ gdcmap3 input.pa3 output.dcm
```

or if you prefer being explicit:

```
$ gdcmap3 -i input.pa3 -o output.dcm
```

Even though **gdcmap3** can overwrite directly on the same file (input.pa3 = output.dcm), it is recommended that user should first convert into a different file to make sure the bug is properly handled by GDCM.

### 11.6 SEE ALSO

**gdcmdump**(1), **gdcmap3**(1), **gdcminfo**(1)

### 11.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 12

# Tool to convert PDF to PDF/DICOM.

### 12.1 SYNOPSIS

```
gdcmpdf [options] file-in file-out
```

### 12.2 DESCRIPTION

The **gdcmpdf** tool convert a PDF file (any PDF version) into an encapsulated PDF/DICOM file. By default it will try to read the PDF meta information stored in the PDF and convert this information to some specific DICOM fields (see below). However it may fails (eg. wrong password on encrypted PDF file) in which case empty value are used.

### 12.3 PARAMETERS

file-in    PDF input filename

file-out   DICOM output filename

### 12.4 OPTIONS

#### 12.4.1 general options

```
-h    --help  
      print this help text and exit  
  
-v    --version  
      print version information and exit  
  
-V    --verbose  
      verbose mode (warning+error).  
  
-W    --warning  
      warning mode, print warning information  
  
-E    --error  
      error mode, print error information  
  
-D    --debug
```

```
debug mode, print debug information
```

## 12.5 Usage Example

```
$ wget http://gdcm.sourceforge.net/gdcm.pdf
$ gdcmpdf gdcm.pdf gdcm.dcm
```

To re-extract the encapsulated PDF file:

```
$ gdcmrw -i gdcm.dcm -t 42,11 -o gdcm.dcm.pdf
$ diff gdcm.pdf gdcm.dcm.pdf
```

## 12.6 PDF Info Mapping

Here is how the PDF info is mapped to DICOM information (typical pdfinfo output):

```
Title:      GDCM Reference Manual
Subject:    Grassroots DICOM API reference
Keywords:   GDCM,DICOM,JPEG,Lossless JPEG,JPEG-LS,J2K,JPEG 2000,RLE
Author:     Mathieu Malaterre and co.
Creator:    LaTeX with hyperref package
Producer:   pdfTeX-1.21a
CreationDate: Tue Apr 28 15:34:26 2009
Tagged:     no
Pages:      1188
Encrypted:  no
Page size:  612 x 792 pts (letter)
File size:  13756841 bytes
Optimized:  yes
PDF version: 1.4
```

Converted to DICOM this leads to:

```
# Dicom-Data-Set
# Used TransferSyntax: Little Endian Explicit
(0008,0005) CS [ISO_IR 100] # 10, 1 SpecificCharacterSet
(0008,0012) DA [20090428] # 8, 1 InstanceCreationDate
(0008,0013) TM [182550.302631] # 14, 1 InstanceCreationTime
(0008,0016) UI =EncapsulatedPDFStorage # 30, 1 SOPClassUID
(0008,0018) UI [1.2.826.0.1.3680043.2.1143.776842935192792959289022034349197114] # 64, 1 SOPInstanceUID
(0008,0020) DA [20090428] # 8, 1 StudyDate
(0008,0023) DA [20090428] # 8, 1 ContentDate
(0008,002a) DT [20090428153437.000000] # 22, 1 AcquisitionDateTime
(0008,0030) TM [182550.302160] # 14, 1 StudyTime
(0008,0033) TM [153426.000000] # 14, 1 ContentTime
(0008,0050) SH (no value available) # 0, 0 AccessionNumber
(0008,0060) CS [OT] # 2, 1 Modality
(0008,0064) CS [WSD] # 4, 1 ConversionType
(0008,0070) LO [LaTeX with hyperref package] # 28, 1 Manufacturer
(0008,0090) PN (no value available) # 0, 0 ReferringPhysiciansName
(0010,0010) PN [Mathieu Malaterre and co.] # 26, 1 PatientsName
(0010,0020) LO (no value available) # 0, 0 PatientID
(0010,0030) DA (no value available) # 0, 0 PatientsBirthDate
(0010,0040) CS (no value available) # 0, 0 PatientsSex
(0018,1020) LO [pdfTeX-1.21a] # 14, 1 SoftwareVersions
(0020,000d) UI [1.2.826.0.1.3680043.2.1143.1868121832223417351654232480755123133] # 64, 1 StudyInstanceUID
(0020,000e) UI [1.2.826.0.1.3680043.2.1143.1330099150825746617507846107663964311] # 64, 1 SeriesInstanceUID
(0020,0010) SH (no value available) # 0, 0 StudyID
(0020,0011) IS [1] # 2, 1 SeriesNumber
```

```

(0020,0013) IS [1] # 2, 1 InstanceNumber
(0028,0301) CS [YES] # 4, 1 BurnedInAnnotation
(0040,a043) SQ (Sequence with explicit length #=0) # 0, 1 ConceptNameCodeSequence
(ffff,e0dd) na (SequenceDelimitationItem for re-encod.) # 0, 0 SequenceDelimitationItem
(0042,0010) ST [GDCM Reference Manual] # 22, 1 DocumentTitle
(0042,0011) OB 25\\50\\44\\46\\2d\\31\\2e\\34\\0a\\25\\e7\\f3\\cf\\d3\\0a\\33\\32\\30\\37\\37\\20\\30... # 137568
(0042,0012) LO [application/pdf] # 16, 1 MIMETimeTypeOfEncapsulatedDocument

```

```

$ stat gdc.m.pdf
  File: `gdc.m.pdf'
  Size: 13756841      Blocks: 26912      IO Block: 4096   regular file
Device: fe01h/65025d Inode: 2675750      Links: 1
Access: (0644/-rw-r--r--)  Uid: ( 1002/mmalaterre)   Gid: ( 1002/mmalaterre)
Access: 2009-04-28 16:05:00.000000000 +0200
Modify: 2009-04-28 15:34:37.000000000 +0200
Change: 2009-04-28 16:05:00.000000000 +0200

```

Explanation for the different Date/Time mappings:

- Study Date/Time, Instance Creation Date/Time are both equal to the current time gdc.mpdf tool was run,
- Acquisition Date Time is set to the Modify Time of the actual PDF file,
- Content Date/Time are set from the actual PDF header info: CreationDate.

## 12.7 SEE ALSO

**gdc.mconv(1), gdc.mraw(1), pdfinfo(1)**

## 12.8 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 13

# Extract Data Element Value Field.

### 13.1 SYNOPSIS

```
gdcmmraw [options] file-in file-out
```

### 13.2 DESCRIPTION

The **gdcmmraw** tool is mostly used for development purpose. It is used to extract a specific binary field from a DICOM DataSet.

### 13.3 PARAMETERS

```
file-in    DICOM input filename
```

```
file-out    output filename
```

### 13.4 OPTIONS

#### 13.4.1 PARAMETERS

```
-i --input      Input filename
-o --output      Output filename
-t --tag        Specify tag to extract value from.
```

#### 13.4.2 OPTIONS

```
-S --split-frags  Split fragments into multiple files.
-p --pattern      Specify trailing file pattern (see split-frags).
-P --pixel-data   Pixel Data trailing 0.
```

#### 13.4.3 general options

```
-h    --help
```







```
-rw-r--r-- 1 mathieu mathieu 81512 2008-08-08 22:10 jpeg03.ljpeg  
-rw-r--r-- 1 mathieu mathieu 81694 2008-08-08 22:10 jpeg02.ljpeg  
-rw-r--r-- 1 mathieu mathieu 81564 2008-08-08 22:10 jpeg01.ljpeg  
-rw-r--r-- 1 mathieu mathieu 79970 2008-08-08 22:10 jpeg00.ljpeg
```

## 13.6 Footnote about JPEG files

It is a common misunderstanding to interchange 'JPEG 8bits lossy' with simply JPEG file. The JPEG specification is much broader than simply the common lossy 8bits file (as found on internet).

You can have:

- JPEG Lossy 8bits
- JPEG Lossy 12bits
- JPEG Lossless 2-16bits

Those are what is defined in ITU-T T.81, ISO/IEC IS 10918-1.

## 13.7 SEE ALSO

**gdcmdump(1)**, **gdcmrw(1)**

## 13.8 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 14

# Scan a directory containing DICOM files.

### 14.1 SYNOPSIS

```
gdcmscanner [options] directory
```

### 14.2 DESCRIPTION

The **gdcmscanner** is a command line tool to quickly extract value from a set of DICOM attribute in a DICOM File-Set.

#### 14.2.1 PARAMETERS

```
-d --dir          DICOM directory
-t --tag %d,%d    DICOM tag(s) to look for
```

#### 14.2.2 OPTIONS

```
-p --print        Print output.
-r --recursive    Recursively descend directory.
```

#### 14.2.3 general options

```
-h  --help
    print this help text and exit

-v  --version
    print version information and exit

-V  --verbose
    verbose mode (warning+error).

-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

## 14.3 Typical usage

## 14.4 Simple usage

In order to display all the value for Patient Name (0010,0010) in the directory name **gdcData**, simply do:

```
$ gdcscanner -t 10,10 -d gdcData -p
```

## 14.5 Complex usage

Because gdcscanner does not support progress, you have to wait until all files are traversed to see any results. This is quite cumbersome, on UNIX this can be worked around with the following trick:

```
$ find gdcData -type d -exec gdcscanner -t 10,10 -d {} -p ';'`
```

So all directory are locally traversed (no child directory are recursively traversed), which means results comes out much faster.

## 14.6 SEE ALSO

**gdcmdump(1)**, **gdcmrw(1)**

## 14.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 15

# Tool to execute a DICOM Query/Retrieve operation

### 15.1 SYNOPSIS

```
gdcmscu [OPTION]...[OPERATION]...HOSTNAME...[PORT]...
```

Execute a DICOM Q/R operation to HOSTNAME, using port PORT (104 when not specified)

### 15.2 DESCRIPTION

The **gdcmscu** command line program is the tool to execute DICOM Query/Retrieve operation. It supports:

- C-ECHO (SCU)
- C-FIND (SCU)
- C-STORE (SCU)
- C-MOVE (SCU/SCP) C-MOVE operation are executed using two different ports (one for the SCU and one for the SCP).

### 15.3 PARAMETERS

### 15.4 OPTIONS

#### 15.4.1 OPTIONS

```
-H --hostname    %s  Hostname.  
-p --port        %d  Port number.  
    --aetitle    %s  Set calling AE Title.  
    --call       %s  Set called AE Title.
```

#### 15.4.2 mode options

```
--echo          C-ECHO (default when none).  
--store         C-STORE.
```

```
--find      C-FIND.
--move      C-MOVE.
```

### 15.4.3 C-STORE options

```
-i --input      %s  DICOM filename
-r --recursive  recursively process (sub-)directories
--store-query %s  Store constructed query in file
```

### 15.4.4 C-FIND/C-MOVE options

```
--patientroot  C-FIND Patient Root Model.
--studyroot    C-FIND Study Root Model.

--patient      C-FIND Query on Patient Info (cannot be used with --studyroot).
--study        C-FIND Query on Study Info.
--series       C-FIND Query on Series Info.
--image        C-FIND Query on Image Info.
--key %d,%d[%s] 0123,4567=VALUE for specifying search criteria (wildcard allowed)
                With --key, leave blank (ie, --key 10,20="" or --key 10,20) to retrieve values
```

### 15.4.5 C-MOVE options

```
-o --output      %s  DICOM filename / directory
--port-scp %d      Port for incoming associations
--key %d,%d[%s]    0123,4567=VALUE for specifying search criteria (wildcard not allowed)
                Note that C-MOVE supports the same queries as C-FIND, but no wildcards are allowed
```

### 15.4.6 general options

```
-h --help
    print this help text and exit

-v --version
    print version information and exit

-V --verbose
    verbose mode (warning+error).

-W --warning
    warning mode, print warning information

-E --error
    error mode, print error information

-D --debug
    debug mode, print debug information

-L --log-file
    specify a filename where to write logs

--queryhelp
    print query help
```

### 15.4.7 environment variable

```
GDCM_ROOT_UID Root UID
```

## 15.5 C-ECHO usage

**gdcm SCU** is a great tool to test if a DICOM server is up. For example to send a C-ECHO to server `dicom.example.com` using port 104, use:

```
$ gdcm SCU dicom.example.com
```

or if you prefer being explicit:

```
$ gdcm SCU --echo dicom.example.com 104
```

Using basic security your DICOM server might require that you set the appropriate called AE-TITLE

```
$ gdcm SCU --echo dicom.example.com 11112 --call SERVSCP
```

If you want to specify your own AE-TITLE (default is GDCMSCU), simply use:

```
$ gdcm SCU --echo dicom.example.com 11112 --call SERVSCP --aetitle MYSCU
```

For example you could test on the DICOM server provided by DICOMObject team:

```
$ gdcm SCU www.dicomserver.co.uk 11112
```

## 15.6 C-STORE usage

C-STORE is the operation that allow sending a DICOM file to a remote DICOM server. For instance to send a file called `myfile.dcm`

```
$ gdcm SCU --store dicom.example.com 104 myfile.dcm
```

or if you prefer being explicit:

```
$ gdcm SCU --store dicom.example.com 104 -i myfile.dcm
```

You can even send multiple files using the same association:

```
$ gdcm SCU --store dicom.example.com 104 myfile1.dcm myfile2.dcm myfile3.dcm ...
```

## 15.7 C-FIND usage

**gdcm SCU** also allow querying a DICOM server. This is the C-FIND operation, for example to find all DICOM Instance where PatientsName match a particular pattern, usage is simply:

```
$ gdcm SCU --find --patient dicom.example.com 11112 --patientroot --key 10,10,"A*"
```

We also support a DCMTK compatible convention:

```
$ gdcm SCU --find --patient dicom.example.com 11112 --patientroot --key 10,10="A*"
```

When an attribute is set without a value it will be part of the output result:

```
$ gdcm SCU --find --patient dicom.example.com 11112 --call MI2B2 --patientroot -k 10,10="A*" -k 10,20
```

## 15.8 C-MOVE usage

C-MOVE is the operation to retrieve a DICOM instance from a remote DICOM server. Most of the time, it is a subsequent operation after a C-FIND query. To retrieve a DICOM instance where PatientID is ABCD1234, simply execute:

```
$ gdcmscu --move --patient --aetitle ACME1 --call ACME_STORE dicom.example.com 5678 --patientroot -k 10,20="ABCD1234"
```

**WARNING** For this operation to work you need information from the DICOM server you are communicating with. Only the DICOM server you are sending a C-MOVE query will be responsible for sending back incoming associations (the actual C-STORE SCP). Therefore you need to make sure that you mapping of (AE-TITLE,PortNumber) is properly set on the DICOM server side as well as the port for incoming association (`--port-scp`).

**gdcmscu** does not currently support external C-STORE association (C-STORE request sent to an external SCP application).

## 15.9 patientroot notes

The flag `--patientroot` is just simply a wrapper around the syntax `--key 8,52=PATIENT` For instance one would write using DCMTK syntax:

```
$ findscu --patient dicom.example.com 11112 --key 8,52=PATIENT --key 10,10="F*"
```

This would become using GDCM syntax:

```
$ gdcmscu --find --patient dicom.example.com 11112 --patientroot --key 10,10="F*"
```

## 15.10 Debugging

This is sometime difficult to investigate why a connection to a remote DICOM server cannot be done. Some recommendations follows:

Always try to do a simple C-ECHO at first. If you cannot get the C-ECHO to work none of the other operations will work. Before trying to a C-MOVE operation, make sure you can execute the C-FIND equivalent query first.

When doing a C-MOVE operation you really need to communicate with the PACS admin as the C-MOVE operation is different from the other lower level operation such as HTTP/GET. When doing a C-MOVE, the server will communicate back using another channel (could be different port) using it's internal database to map an AE-TITLE back to the destination IP.

Indeed the C-MOVE operation by design does not always use your incoming IP address to send back the resulting dataset. Instead it uses a mapping of AE-TITLE to IP address to send back any results. So pay particular attention to the spelling of your AE-TITLE and your incoming port (which may be different from the port to connect to the server).

## 15.11 Port Warning

Watch out that port ranging [1-1024] are reserved for admin and not easily accessible unless granted special privileges. Therefore the default 104 DICOM port might not be accessible to all your users.



## 15.12 C-STORE Warnings

When constructing a C-STORE operation, `gdcm SCU` will always use the Media Storage SOP Class UID as found in the file to be sent. For encapsulated DICOM file (eg. RLE Lossless) the receiving SCP server might not support this compression and will legitimately refuse the C-STORE operation. In this case users have to manually convert to a non-compressed form this particular file:

```
$ gdcmconv --raw compressed.dcm non_compressed.dcm
```

## 15.13 C-MOVE Warnings

At the moment `gdcm SCU` only supports non-compressed transfer syntax. It will always request DataSet using Implicit VR Little Endian Transfer Syntax during a C-MOVE operation (both incoming and outgoing associations). This make `gdcm SCU -move` equivalent to DCMTK `movescu` syntax:

```
$ movescu -xi +xi ...
```

## 15.14 C-FIND IMAGE level (Composite Object Instance)

One should pay attention that `gdcm SCU -find` and `findscu` are not completely equivalent. Using `gdcm SCU -find`, all Unique Keys will be added automatically. One can therefore execute something like this:

```
$ gdcm SCU --find --patientroot --image --key 8,18=1.2.3.4.5.6 dicom.example.com 11112
```

instead of the more explicit form

```
$ gdcm SCU --find --patientroot --image --key 8,18=1.2.3.4.5.6 dicom.example.com 11112 --key 10,20 --key 20,d --key
```

This would also be equivalent to:

```
$ findscu --patient --key 8,52=IMAGE --key 8,18=1.2.3.4.5.6 dicom.example.com 11112 --key 10,20 --key 20,d --key
```

## 15.15 Storing the Query

It is also possible to store the query:

```
gdcm SCU --find --patient --patientroot dicom.example.com 11112 --key 10,20="*" --key 10,10 --store-query query.dcm
```

One can then check the DataSet values send for the query:

```
$ gdcmdump query.dcm
# Dicom-File-Format

# Dicom-Meta-Information-Header
# Used TransferSyntax:

# Dicom-Data-Set
# Used TransferSyntax: 1.2.840.10008.1.2
(0008,0005) ?? (CS) [ISO_IR 192] # 10,1-n Specific Character Set
(0008,0052) ?? (CS) [PATIENT ] # 8,1 Query/Retrieve Level
(0010,0010) ?? (PN) (no value) # 0,1 Patient's Name
(0010,0020) ?? (LO) [* ] # 2,1 Patient ID
```

The Specific Character Set was set to "ISO\_IR 192" as the locale encoding of the system was found automatically by gdcmscu to be UTF-8.

This means that the following command line will properly setup the Query with the appropriate Charset to be executed correctly:

```
$ gdcmscu --find --patient --patientroot dicom.example.com 11112 --key 10,10="*Jérôme"
```

The query is always executed on the server side (SCP), some implementations does not support string matching with different Character Set.

## 15.16 DICOM Public Servers

An up to date list of DICOM Public Servers can be found at:

<http://www.dclunie.com/medical-image-faq/html/part8.html#DICOMPublicServers>

## 15.17 SEE ALSO

**gdcconv(1)**

## 15.18 COPYRIGHT

Copyright Insight Software Consortium

## Chapter 16

# Concatenate/Extract DICOM files.

### 16.1 SYNOPSIS

```
gdcmtar [options] file-in file-out
```

### 16.2 DESCRIPTION

The **gdcmtar** is a command line tool used to tar/untar multi-frames images (including SIEMENS MOSAIC file)

### 16.3 PARAMETERS

file-in    DICOM input filename

file-out   DICOM output filename

### 16.4 OPTIONS

#### 16.4.1 OPTIONS

```
--enhance      enhance (default)
-U --unenhance  unenhance
-M --mosaic     Split SIEMENS Mosaic image into multiple frames.
-p --pattern    Specify trailing file pattern.
--root-uid      Root UID.
```

#### 16.4.2 general options

```
-h --help      print this help text and exit
-v --version    print version information and exit
-V --verbose    verbose mode (warning+error).
```

```
-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

### 16.4.3 environment variable

GDCM\_ROOT\_UID Root UID

## 16.5 Typical usage

### 16.5.1 SIEMENS Mosaic

```
$ gdcminfo MR-sonata-3D-as-Tile.dcm
```

```
MediaStorage is 1.2.840.10008.5.1.4.1.1.4 [MR Image Storage]
TransferSyntax is 1.2.840.10008.1.2.1 [Explicit VR Little Endian]
NumberOfDimensions: 2
Dimensions: (384,384,1)
\&...
```

```
$ gdcmtar --mosaic -i MR-sonata-3D-as-Tile.dcm -o mosaic --pattern %03d.dcm
```

Will output:

```
-rw-r--r-- 1 mathieu mathieu 72882 2009-08-10 11:14 mosaic000.dcm
-rw-r--r-- 1 mathieu mathieu 72886 2009-08-10 11:14 mosaic001.dcm
-rw-r--r-- 1 mathieu mathieu 72886 2009-08-10 11:14 mosaic002.dcm
-rw-r--r-- 1 mathieu mathieu 72886 2009-08-10 11:14 mosaic003.dcm
-rw-r--r-- 1 mathieu mathieu 72886 2009-08-10 11:14 mosaic004.dcm
-rw-r--r-- 1 mathieu mathieu 72886 2009-08-10 11:14 mosaic005.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic006.dcm
-rw-r--r-- 1 mathieu mathieu 72882 2009-08-10 11:14 mosaic007.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic008.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic009.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic010.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic011.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic012.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic013.dcm
-rw-r--r-- 1 mathieu mathieu 72882 2009-08-10 11:14 mosaic014.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic015.dcm
-rw-r--r-- 1 mathieu mathieu 72882 2009-08-10 11:14 mosaic016.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic017.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic018.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic019.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic020.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic021.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic022.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic023.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic024.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic025.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic026.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic027.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic028.dcm
-rw-r--r-- 1 mathieu mathieu 72884 2009-08-10 11:14 mosaic029.dcm
-rw-r--r-- 1 mathieu mathieu 72882 2009-08-10 11:14 mosaic030.dcm
```

```
$ gdcminfo mosaic000.dcm
```

```
MediaStorage is 1.2.840.10008.5.1.4.1.1.4 [MR Image Storage]  
TransferSyntax is 1.2.840.10008.1.2.1 [Explicit VR Little Endian]  
NumberOfDimensions: 2  
Dimensions: (64,64,1)  
\&...
```

## 16.6 SEE ALSO

**gdcmdump(1)**, **gdcmrw(1)**, **gdcminfo(1)**

## 16.7 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 17

# Simple DICOM viewer.

### 17.1 SYNOPSIS

```
gdcviewer [options] file-in
```

### 17.2 DESCRIPTION

The **gdcviewer** is a simple tool that show how to use [vtkGDCMImageReader](#). The class that use gdc to make a layer to VTK. **gdcviewer** is basically only just a wrapper around VTK/GDCM.

This tool is meant for testing integration of GDCM in VTK. You should see it as a demo tool. It does compile with VTK ranging from 4.2 to 5.10, but only with VTK 5.2 (or above) can play with the widgets (as described below).

### 17.3 PARAMETERS

```
file-in    DICOM input filename
```

### 17.4 OPTIONS

#### 17.4.1 OPTIONS

<code>--force-rescale</code>	force rescale (advanced users)
<code>--force-spacing</code>	force spacing (advanced users)
<code>-r --recursive</code>	Recursively descend directory

#### 17.4.2 general options

<code>-h</code>	<code>--help</code>	print this help text and exit
<code>-v</code>	<code>--version</code>	print version information and exit
<code>-V</code>	<code>--verbose</code>	verbose mode (warning+error).

```
-W  --warning
    warning mode, print warning information

-E  --error
    error mode, print error information

-D  --debug
    debug mode, print debug information
```

## 17.5 Typical usage

## 17.6 Simple usage

For now `gdcmviewer` should be started from a command line prompt. The next argument should be the name of the DICOM file you wish to read. For instance:

```
$ gdcmviewer -V 012345.002.050.dcm
```

`gdcmviewer` will try to read your file, and then print the `vtk` information associated with this file. Basically what kind of image you are looking at.

- `ScalarType` is the DICOM Real World Value type
- `Dimensions` is the dimension of the image
- `Spacing` is the spacing of the image
- `NumberOfScalarComponents` should be 1 for grayscale & `PALETTE COLOR` and 3 for `RGB`, `YBR` data.

## 17.7 Wiki Link

The wiki page, with color pictures can be found at: <http://gdcm.sourceforge.net/wiki/index.php/Gdcmviewer>

## 17.8 SEE ALSO

`gdcmdump(1)`, `gdcm2vtk(1)`

## 17.9 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre



## Chapter 18

**provides a tool to convert a DICOM file into a XML info set and vice-versa.**

### 18.1 SYNOPSIS

```
gdcxml [options] file-in[DICOM or XML] file-out[XML or DICOM]
```

### 18.2 DESCRIPTION

The **gdcxml** command line program converts a DICOM file (DataSet) into an XML file (according to the Native DICOM Model) or vice-versa. For those familiar with DCMTK, this provides binary capabilities (i.e. functionality of both dcm2xml and xml2dcm).

The XML info set which is from the DICOM file gdcXMLPrintet Class. This is in strict compliance with the Native DICOM Model as given in Supp 118.

### 18.3 PARAMETERS

```
file-in    DICOM or XML input filename ( cannot be absent)
```

```
file-out    output filename (can be absent)
```

### 18.4 OPTIONS

#### 18.4.1 PARAMETERS

```
-i --input      DICOM filename
-o --output      DICOM filename
```

#### 18.4.2 Options for DICOM to XML:

```
-B --loadBulkData  Loads bulk data into a binary file named "UUID" (by default UUID are written).
```

### 18.4.3 Options for XML to DICOM:

```
-B --loadBulkData  Loads bulk data from a binary file named as the "UUID" in XML file (by default UUID are writ  
-T --TransferSyntax Loads transfer syntax from file (default is LittleEndianImplicit)
```

### 18.4.4 general options

```
-h  --help  
    print this help text and exit  
  
-v  --version  
    print version information and exit  
  
-V  --verbose  
    verbose mode (warning+error).  
  
-W  --warning  
    warning mode, print warning information  
  
-E  --error  
    error mode, print error information  
  
-D  --debug  
    debug mode, print debug information
```

## 18.5 SEE ALSO

**gdcmdump(1), gdcconv(1)**

## 18.6 COPYRIGHT

Copyright (c) 2006-2011 Mathieu Malaterre

## Chapter 19

# Todo List

### Class `gdcm::CSAHeader`

MrEvaProtocol in 29,1020 contains ^M that would be nice to get rid of on UNIX system...

### Class `gdcm::network::ApplicationContext`

Looks like Application Context can only be 64 bytes at max (see Figure 9-1 / PS 3.8 - 2009 )

### Class `gdcm::Overlay`

Is there actually any way to recognize an overlay ? On images with multiple overlay I do not see any way to differentiate them (other than the group tag).

### Class `gdcm::SequenceOfFragments`

I do not enforce that Sequence of Fragments ends with a SQ end del

### Class `gdcm::TransferSyntax`

: The implementation is completely retarded -> see `gdcm::UIDs` for a replacement We need: `IsSupported` We need preprocess of raw/xml file We need `GetFullName()`

### Member `gdcm::UIDGenerator::IsValid (const char *uid)`

: Move that in `DataStructureAndEncoding` (see `FileMetaInformation::CheckFileMetaInformation`)



## Chapter 20

# Deprecated List

**Member [gdcm::CompositeNetworkFunctions::ConstructQuery](#)** (ERootType inRootType, EQueryLevel inQueryLevel, const KeyValuePairArrayType &keys, bool inMove=false)

**Member [gdcm::FileSet::AddFile](#)** (File const &)

. Does nothing

**Member [gdcm::TransferSyntax::GetSwapCode](#)** () const

Return the [SwapCode](#) associated with the Transfer Syntax. Be careful with the special GE private syntax the [DataSet](#) is written in little endian but the Pixel Data is in Big Endian.



## Chapter 21

# Bug List

### Class `gdcm::DICOMDIRGenerator`

: There is a current limitation of not handling Referenced SOP Class UID / Referenced SOP Instance UID simply because the `Scanner` does not allow us See PS 3.11 / [Table D.3-2 STD-GEN Additional DICOMDIR Keys](#)

### Class `gdcm::IPPSorter`

There are currently a couple of bugs in this implementation:

- Gantry Tilt is not considered





## Chapter 22

# Namespace Index

### 22.1 Namespace List

Here is a list of all namespaces with brief descriptions:

<a href="#">gdc</a>	109
<a href="#">gdc::network</a>	131
<a href="#">gdc::SegmentHelper</a>	136
<a href="#">gdc::terminal</a>	
Class for Terminal Allow one to print in color in a shell	136



## Chapter 23

# Hierarchical Index

### 23.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

gdcn::network::AbstractSyntax	150
gdcn::network::ApplicationContext	160
gdcn::ApplicationEntity	161
gdcn::network::ARTIMTimer	166
gdcn::ASN1	167
gdcn::network::AsynchronousOperationsWindowSub	168
gdcn::Attribute< Group, Element, TVR, TVM >	169
gdcn::Attribute< Group, Element, TVR, VM::VM1 >	176
gdcn::Attribute< Group, Element, TVR, VM::VM1_n >	183
gdcn::Attribute< Group, Element, TVR, VM::VM1_3 >	180
gdcn::Attribute< Group, Element, TVR, VM::VM1_8 >	181
gdcn::Attribute< Group, Element, TVR, VM::VM2_n >	188
gdcn::Attribute< Group, Element, TVR, VM::VM2_2n >	187
gdcn::Attribute< Group, Element, TVR, VM::VM3_n >	191
gdcn::Attribute< Group, Element, TVR, VM::VM3_3n >	190
gdcn::Base64	195
gdcn::network::BaseCompositeMessage	196
gdcn::network::CEchoRQ	233
gdcn::network::CEchoRSP	234
gdcn::network::CFindCancelRQ	236
gdcn::network::CFindRQ	237
gdcn::network::CFindRSP	238
gdcn::network::CMoveCancelRq	239
gdcn::network::CMoveRQ	241
gdcn::network::CMoveRSP	242
gdcn::network::CStoreRQ	277
gdcn::network::CStoreRSP	278
gdcn::network::BasePDU	198
gdcn::network::AAabortPDU	139
gdcn::network::AAssociateACPDU	141
gdcn::network::AAssociateRJPDU	144
gdcn::network::AAssociateRQPDU	145
gdcn::network::AReleaseRPPDU	163

gdcmm::network::AReleaseRQPDU . . . . .	164
gdcmm::network::PDataTFPDU . . . . .	556
std::basic_string< Char >	
std::string	
gdcmm::String< TDelimiter, TMaxLength, TPadChar > . . . . .	707
gdcmm::SegmentHelper::BasicCodedEntry . . . . .	204
gdcmm::BitmapToBitmapFilter . . . . .	217
gdcmm::PixmapToPixmapFilter . . . . .	583
gdcmm::ImageToImageFilter . . . . .	449
gdcmm::ImageApplyLookupTable . . . . .	417
gdcmm::ImageChangePhotometricInterpretation . . . . .	419
gdcmm::ImageChangePlanarConfiguration . . . . .	422
gdcmm::ImageChangeTransferSyntax . . . . .	425
gdcmm::ImageFragmentSplitter . . . . .	436
gdcmm::ByteBuffer . . . . .	222
gdcmm::ByteSwap< T > . . . . .	223
gdcmm::ByteSwapFilter . . . . .	224
gdcmm::network::CFind . . . . .	235
gdcmm::Coder . . . . .	244
gdcmm::Codec . . . . .	243
gdcmm::AudioCodec . . . . .	193
gdcmm::ImageCodec . . . . .	429
gdcmm::DeltaEncodingCodec . . . . .	309
gdcmm::JPEG2000Codec . . . . .	476
gdcmm::JPEGCodec . . . . .	482
gdcmm::JPEG12Codec . . . . .	471
gdcmm::JPEG16Codec . . . . .	474
gdcmm::JPEG8Codec . . . . .	479
gdcmm::JPEGLSCCodec . . . . .	486
gdcmm::KAKADUCCodec . . . . .	491
gdcmm::PGXCodec . . . . .	567
gdcmm::PNMCodec . . . . .	588
gdcmm::PVRGCodec . . . . .	610
gdcmm::RAWCodec . . . . .	624
gdcmm::RLECodec . . . . .	637
gdcmm::PDFCodec . . . . .	562
gdcmm::CodeString . . . . .	245
gdcmm::network::CompositeMessageFactory . . . . .	252
gdcmm::CompositeNetworkFunctions . . . . .	253
gdcmm::ConstCharWrapper . . . . .	258
gdcmm::CryptoFactory . . . . .	260
gdcmm::CAPICryptoFactory . . . . .	229
gdcmm::OpenSSLCryptoFactory . . . . .	537
gdcmm::OpenSSL7CryptoFactory . . . . .	540
gdcmm::CryptographicMessageSyntax . . . . .	262
gdcmm::CAPICryptographicMessageSyntax . . . . .	230
gdcmm::OpenSSLCryptographicMessageSyntax . . . . .	538
gdcmm::OpenSSL7CryptographicMessageSyntax . . . . .	542
gdcmm::CSAElement . . . . .	264
gdcmm::CSAHeader . . . . .	268
gdcmm::CSAHeaderDict . . . . .	272
gdcmm::CSAHeaderDictEntry . . . . .	274

gdcm::DataElement . . . . .	282
gdcm::CP246ExplicitDataElement . . . . .	258
gdcm::ExplicitDataElement . . . . .	362
gdcm::ExplicitImplicitDataElement . . . . .	363
gdcm::Fragment . . . . .	401
gdcm::BasicOffsetTable . . . . .	207
gdcm::ImplicitDataElement . . . . .	456
gdcm::Item . . . . .	467
gdcm::UNExplicitDataElement . . . . .	830
gdcm::UNExplicitImplicitDataElement . . . . .	831
gdcm::VR16ExplicitDataElement . . . . .	855
gdcm::DataSet . . . . .	294
gdcm::CommandDataSet . . . . .	250
gdcm::FileMetaInformation . . . . .	379
gdcm::DataSetHelper . . . . .	304
gdcm::Decoder . . . . .	304
gdcm::Codec . . . . .	243
gdcm::DefinedTerms . . . . .	306
gdcm::Defs . . . . .	306
gdcm::DICOMDIR . . . . .	310
gdcm::DICOMDIRGenerator . . . . .	311
gdcm::Dict . . . . .	313
gdcm::DictConverter . . . . .	316
gdcm::DictEntry . . . . .	318
gdcm::Dicts . . . . .	322
gdcm::network::DIMSE . . . . .	325
gdcm::DirectionCosines . . . . .	326
gdcm::Directory . . . . .	328
gdcm::DirectoryHelper . . . . .	331
gdcm::DummyValueGenerator . . . . .	332
gdcm::Element< TVR, TVM > . . . . .	335
gdcm::Element< TVR, VM::VM1_n > . . . . .	339
gdcm::Element< TVR, VM::VM1_2 > . . . . .	338
gdcm::Element< TVR, VM::VM2_n > . . . . .	343
gdcm::Element< TVR, VM::VM2_2n > . . . . .	341
gdcm::Element< TVR, VM::VM3_n > . . . . .	346
gdcm::Element< TVR, VM::VM3_3n > . . . . .	344
gdcm::Element< VR::AS, VM::VM5 > . . . . .	347
gdcm::Element< VR::OB, VM::VM1_n > . . . . .	335
gdcm::Element< VR::OB, VM::VM1 > . . . . .	348
gdcm::Element< VR::OW, VM::VM1_n > . . . . .	335
gdcm::Element< VR::OW, VM::VM1 > . . . . .	349
gdcm::ElementDisableCombinations< TVR, TVM > . . . . .	351
gdcm::ElementDisableCombinations< VR::OB, VM::VM1_n > . . . . .	352
gdcm::ElementDisableCombinations< VR::OW, VM::VM1_n > . . . . .	352
gdcm::EncapsulatedDocument . . . . .	352
gdcm::EncodingImplementation< T > . . . . .	353
gdcm::EncodingImplementation< VR::VRASCII > . . . . .	353
gdcm::EncodingImplementation< VR::VRBINARY > . . . . .	354
gdcm::EnumeratedValues . . . . .	356
gdcm::Event . . . . .	357
gdcm::AnyEvent . . . . .	158

gdcmm::AbortEvent . . . . .	149
gdcmm::AnonymizeEvent . . . . .	151
gdcmm::DataEvent . . . . .	292
gdcmm::DataSetEvent . . . . .	301
gdcmm::EndEvent . . . . .	355
gdcmm::ExitEvent . . . . .	361
gdcmm::FileNameEvent . . . . .	386
gdcmm::InitializeEvent . . . . .	458
gdcmm::IterationEvent . . . . .	470
gdcmm::ModifiedEvent . . . . .	519
gdcmm::ProgressEvent . . . . .	608
gdcmm::StartEvent . . . . .	697
gdcmm::UserEvent . . . . .	836
gdcmm::NoEvent . . . . .	533
std::exception	
gdcmm::CSAHeaderDictException . . . . .	276
gdcmm::DataElementException . . . . .	291
gdcmm::Exception . . . . .	359
gdcmm::ParseException . . . . .	552
gdcmm::Fiducials . . . . .	365
gdcmm::FileDerivation . . . . .	374
gdcmm::FileExplicitFilter . . . . .	377
gdcmm::Filename . . . . .	384
gdcmm::FilenameGenerator . . . . .	388
gdcmm::FileSet . . . . .	391
gdcmm::Global . . . . .	404
gdcmm::GroupDict . . . . .	406
gdcmm::IconImageFilter . . . . .	408
gdcmm::IconImageGenerator . . . . .	410
gdcmm::ignore_char . . . . .	412
gdcmm::ImageConverter . . . . .	435
gdcmm::ImageHelper . . . . .	439
gdcmm::network::ImplementationClassUIDSub . . . . .	454
gdcmm::network::ImplementationUIDSub . . . . .	454
gdcmm::network::ImplementationVersionNameSub . . . . .	455
gdcmm::IOD . . . . .	459
gdcmm::IODEntry . . . . .	460
gdcmm::IODs . . . . .	462
gdcmm::JSON . . . . .	490
gdcmm::Scanner::ltstr . . . . .	500
gdcmm::Macro . . . . .	501
gdcmm::Macros . . . . .	502
gdcmm::network::MaximumLengthSub . . . . .	503
gdcmm::MD5 . . . . .	504
gdcmm::MediaStorage . . . . .	505
gdcmm::Module . . . . .	521
gdcmm::ModuleEntry . . . . .	522
gdcmm::NestedModuleEntries . . . . .	531
gdcmm::Modules . . . . .	525
gdcmm::Object . . . . .	534
gdcmm::BaseRootQuery . . . . .	200
gdcmm::FindPatientRootQuery . . . . .	397
gdcmm::FindStudyRootQuery . . . . .	399

gdcmm::MovePatientRootQuery . . . . .	527
gdcmm::MoveStudyRootQuery . . . . .	529
gdcmm::Bitmap . . . . .	209
gdcmm::Pixmap . . . . .	577
gdcmm::Image . . . . .	413
gdcmm::Curve . . . . .	280
gdcmm::File . . . . .	366
gdcmm::FileWithName . . . . .	396
gdcmm::LookupTable . . . . .	496
gdcmm::SegmentedPaletteColorLookupTable . . . . .	652
gdcmm::MeshPrimitive . . . . .	516
gdcmm::Overlay . . . . .	546
gdcmm::Segment . . . . .	648
gdcmm::Subject . . . . .	713
gdcmm::Anonymizer . . . . .	154
gdcmm::Command . . . . .	248
gdcmm::MemberCommand< T > . . . . .	512
gdcmm::SimpleMemberCommand< T > . . . . .	680
gdcmm::FileAnonymizer . . . . .	369
gdcmm::FileChangeTransferSyntax . . . . .	372
gdcmm::FileStreamer . . . . .	392
gdcmm::network::ULConnectionManager . . . . .	824
gdcmm::Scanner . . . . .	642
gdcmm::ServiceClassUser . . . . .	674
gdcmm::Surface . . . . .	716
gdcmm::Value . . . . .	840
gdcmm::ByteValue . . . . .	224
gdcmm::SequenceOfFragments . . . . .	659
gdcmm::SequenceOfItems . . . . .	664
gdcmm::Orientation . . . . .	544
gdcmm::Parser . . . . .	554
gdcmm::Patient . . . . .	556
gdcmm::PDBelement . . . . .	559
gdcmm::PDBHeader . . . . .	561
gdcmm::network::PDUFactory . . . . .	564
gdcmm::PersonName . . . . .	566
gdcmm::PhotometricInterpretation . . . . .	569
gdcmm::PixelFormat . . . . .	571
gdcmm::Preamble . . . . .	591
gdcmm::PresentationContext . . . . .	592
gdcmm::network::PresentationContextAC . . . . .	594
gdcmm::PresentationContextGenerator . . . . .	595
gdcmm::network::PresentationContextRQ . . . . .	597
gdcmm::network::PresentationDataValue . . . . .	599
gdcmm::Printer . . . . .	601
gdcmm::DictPrinter . . . . .	320
gdcmm::Dumper . . . . .	333
gdcmm::PrivateDict . . . . .	604
gdcmm::PythonFilter . . . . .	612
gdcmm::QueryBase . . . . .	613
gdcmm::QueryImage . . . . .	616
gdcmm::QueryPatient . . . . .	618
gdcmm::QuerySeries . . . . .	620

gdcmm::QueryStudy . . . . .	622
gdcmm::QueryFactory . . . . .	615
gdcmm::Reader . . . . .	627
gdcmm::PixmapReader . . . . .	579
gdcmm::ImageReader . . . . .	443
gdcmm::ImageRegionReader . . . . .	446
gdcmm::SegmentReader . . . . .	654
gdcmm::SurfaceReader . . . . .	725
gdcmm::Region . . . . .	632
gdcmm::BoxRegion . . . . .	219
gdcmm::Rescaler . . . . .	634
gdcmm::network::RoleSelectionSub . . . . .	640
gdcmm::SerieHelper::Rule . . . . .	641
gdcmm::SerieHelper . . . . .	670
gdcmm::Series . . . . .	672
gdcmm::network::ServiceClassApplicationInformation . . . . .	673
gdcmm::SHA1 . . . . .	679
gdcmm::SimpleSubjectWatcher . . . . .	684
gdcmm::SmartPointer< ObjectType > . . . . .	685
gdcmm::SmartPointer< gdcmm::Bitmap > . . . . .	685
gdcmm::SmartPointer< gdcmm::File > . . . . .	685
gdcmm::SmartPointer< gdcmm::Image > . . . . .	685
gdcmm::SmartPointer< gdcmm::MemberCommand > . . . . .	685
gdcmm::SmartPointer< gdcmm::MeshPrimitive > . . . . .	685
gdcmm::SmartPointer< gdcmm::Pixmap > . . . . .	685
gdcmm::SmartPointer< gdcmm::SimpleMemberCommand > . . . . .	685
gdcmm::SmartPointer< gdcmm::Subject > . . . . .	685
gdcmm::SmartPointer< LookupTable > . . . . .	685
gdcmm::SmartPointer< Segment > . . . . .	685
gdcmm::SmartPointer< Surface > . . . . .	685
gdcmm::SmartPointer< Value > . . . . .	685
gdcmm::network::SOPClassExtendedNegociationSub . . . . .	688
gdcmm::SOPClassUIDToIOD . . . . .	689
gdcmm::Sorter . . . . .	690
gdcmm::IPPSorter . . . . .	464
gdcmm::Spacing . . . . .	694
gdcmm::Spectroscopy . . . . .	696
gdcmm::SplitMosaicFilter . . . . .	696
gdcmm::static_assert_test< x > . . . . .	699
gdcmm::STATIC_ASSERTION_FAILURE< x > . . . . .	699
gdcmm::STATIC_ASSERTION_FAILURE< true > . . . . .	699
gdcmm::StreamImageReader . . . . .	699
gdcmm::StreamImageWriter . . . . .	702
String< '\', 64 > . . . . .	
gdcmm::LO . . . . .	493
gdcmm::StringFilter . . . . .	710
gdcmm::Study . . . . .	713
gdcmm::SurfaceHelper . . . . .	722
gdcmm::SwapCode . . . . .	729
gdcmm::SwapperDoOp . . . . .	731
gdcmm::SwapperNoOp . . . . .	732
gdcmm::System . . . . .	732
gdcmm::Table . . . . .	736



gdcm::TableEntry . . . . .	737
gdcm::TableReader . . . . .	738
gdcm::XMLDictReader . . . . .	921
gdcm::XMLPrivateDictReader . . . . .	926
gdcm::network::TableRow . . . . .	740
gdcm::Tag . . . . .	741
gdcm::PrivateTag . . . . .	606
gdcm::TagPath . . . . .	748
gdcm::Testing . . . . .	749
gdcm::Trace . . . . .	754
gdcm::TransferSyntax . . . . .	757
gdcm::network::TransferSyntaxSub . . . . .	761
gdcm::network::Transition . . . . .	762
gdcm::Type . . . . .	764
gdcm::UI . . . . .	765
gdcm::UIDGenerator . . . . .	766
gdcm::UIDs . . . . .	768
gdcm::network::ULAction . . . . .	786
gdcm::network::ULActionAA1 . . . . .	789
gdcm::network::ULActionAA2 . . . . .	790
gdcm::network::ULActionAA3 . . . . .	791
gdcm::network::ULActionAA4 . . . . .	792
gdcm::network::ULActionAA5 . . . . .	793
gdcm::network::ULActionAA6 . . . . .	794
gdcm::network::ULActionAA7 . . . . .	795
gdcm::network::ULActionAA8 . . . . .	796
gdcm::network::ULActionAE1 . . . . .	797
gdcm::network::ULActionAE2 . . . . .	798
gdcm::network::ULActionAE3 . . . . .	799
gdcm::network::ULActionAE4 . . . . .	800
gdcm::network::ULActionAE5 . . . . .	801
gdcm::network::ULActionAE6 . . . . .	802
gdcm::network::ULActionAE7 . . . . .	803
gdcm::network::ULActionAE8 . . . . .	804
gdcm::network::ULActionAR1 . . . . .	805
gdcm::network::ULActionAR10 . . . . .	806
gdcm::network::ULActionAR2 . . . . .	807
gdcm::network::ULActionAR3 . . . . .	808
gdcm::network::ULActionAR4 . . . . .	809
gdcm::network::ULActionAR5 . . . . .	810
gdcm::network::ULActionAR6 . . . . .	811
gdcm::network::ULActionAR7 . . . . .	812
gdcm::network::ULActionAR8 . . . . .	813
gdcm::network::ULActionAR9 . . . . .	814
gdcm::network::ULActionDT1 . . . . .	815
gdcm::network::ULActionDT2 . . . . .	816
gdcm::network::ULConnection . . . . .	818
gdcm::network::ULConnectionCallback . . . . .	821
gdcm::network::ULBasicCallback . . . . .	817
gdcm::network::ULWritingCallback . . . . .	828
gdcm::network::ULConnectionInfo . . . . .	822
gdcm::network::ULEvent . . . . .	826
gdcm::network::ULTransitionTable . . . . .	827

gdcm::Unpacker12Bits	833
gdcm::Usage	834
gdcm::network::UserInformation	837
gdcm::UUIDGenerator	838
gdcm::Validate	839
gdcm::ValueIO< TDE, TSwap, TType >	842
gdcm::Version	843
gdcm::VL	844
gdcm::VM	846
gdcm::VMToLength< T >	850
gdcm::VR	850
gdcm::VRToEncoding< T >	857
gdcm::VRToType< T >	857
gdcm::VRToType< TVR >	857
gdcm::VRVLSIZE< T >	858
gdcm::VRVLSIZE< 0 >	858
gdcm::VRVLSIZE< 1 >	858
vtkImageAlgorithm	
vtkImagePlanarComponentsToComponents	904
vtkImageMapToColors	
vtkImageMapToWindowLevelColors2	902
vtkImageWriter	
vtkGDCMImageWriter	871
vtkLookupTable	
vtkLookupTable16	909
vtkMedicalImageProperties	
vtkGDCMMedicalImageProperties	875
vtkMedicalImageReader2	
vtkGDCMImageReader	859
vtkGDCMThreadedImageReader	886
vtkGDCMImageReader2	865
vtkObject	
vtkGDCMTesting	883
vtkImageColorViewer	892
vtkRTStructSetProperties	911
vtkPolyDataAlgorithm	
vtkGDCMPolyDataReader	877
vtkPolyDataWriter	
vtkGDCMPolyDataWriter	880
vtkThreadedImageAlgorithm	
vtkGDCMThreadedImageReader2	888
vtkImageMapToColors16	899
vtkImageRGBToYBR	906
vtkImageYBRToRGB	908
gdcm::Waveform	916
gdcm::Writer	916
gdcm::PixmapWriter	585
gdcm::ImageWriter	451
gdcm::SegmentWriter	656
gdcm::SurfaceWriter	727
gdcm::XMLPrinter	923

## Chapter 24

# Class Index

### 24.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

<a href="#">gdcn::network::AAabortPDU</a>	
<a href="#">AAabortPDU Table</a>	9-26 A-ABORT PDU FIELDS . . . . . 139
<a href="#">gdcn::network::AAssociateACPDU</a>	
<a href="#">AAssociateACPDU Table</a>	9-17 ASSOCIATE-AC PDU fields . . . . . 141
<a href="#">gdcn::network::AAssociateRJPDU</a>	
<a href="#">AAssociateRJPDU Table</a>	9-21 ASSOCIATE-RJ PDU FIELDS . . . . . 144
<a href="#">gdcn::network::AAssociateRQPDU</a>	
<a href="#">AAssociateRQPDU Table</a>	9-11 ASSOCIATE-RQ PDU fields . . . . . 145
<a href="#">gdcn::AbortEvent</a>	. . . . . 149
<a href="#">gdcn::network::AbstractSyntax</a>	
<a href="#">AbstractSyntax Table</a>	9-14 ABSTRACT SYNTAX SUB-ITEM FIELDS . . . . . 150
<a href="#">gdcn::AnonymizeEvent</a>	
<a href="#">AnonymizeEvent</a>	Special type of event triggered during the Anonymization process . . . . . 151
<a href="#">gdcn::Anonymizer</a>	
<a href="#">Anonymizer</a>	This class is a multi purpose anonymizer. It can work in 2 mode: . . . . . 154
<a href="#">gdcn::AnyEvent</a>	. . . . . 158
<a href="#">gdcn::network::ApplicationContext</a>	
<a href="#">ApplicationContext Table</a>	9-12 APPLICATION CONTEXT ITEM FIELDS . . . . . 160
<a href="#">gdcn::ApplicationEntity</a>	
<a href="#">ApplicationEntity</a>	. . . . . 161
<a href="#">gdcn::network::AReleaseRPPDU</a>	
<a href="#">AReleaseRPPDU Table</a>	9-25 A-RELEASE-RP PDU fields . . . . . 163
<a href="#">gdcn::network::AReleaseRQPDU</a>	
<a href="#">AReleaseRQPDU Table</a>	9-24 A-RELEASE-RQ PDU FIELDS . . . . . 164
<a href="#">gdcn::network::ARTIMTimer</a>	
<a href="#">ARTIMTimer</a>	This file contains the code for the ARTIM timer . . . . . 166
<a href="#">gdcn::ASN1</a>	
Class for <a href="#">ASN1</a>	. . . . . 167
<a href="#">gdcn::network::AsynchronousOperationsWindowSub</a>	
<a href="#">AsynchronousOperationsWindowSub</a>	PS 3.7 <a href="#">Table D.3-7</a> ASYNCHRONOUS OPERATIONS WIN↔
DOW SUB-ITEM FIELDS (A-ASSOCIATE-RQ)	. . . . . 168
<a href="#">gdcn::Attribute&lt; Group, Element, TVR, TVM &gt;</a>	
<a href="#">Attribute</a> class	This class use template metaprograming tricks to let the user know when the template instantiation does not match the public dictionary . . . . . 169

gdcm::Attribute< Group, Element, TVR, VM::VM1 >	176
gdcm::Attribute< Group, Element, TVR, VM::VM1_3 >	180
gdcm::Attribute< Group, Element, TVR, VM::VM1_8 >	181
gdcm::Attribute< Group, Element, TVR, VM::VM1_n >	183
gdcm::Attribute< Group, Element, TVR, VM::VM2_2n >	187
gdcm::Attribute< Group, Element, TVR, VM::VM2_n >	188
gdcm::Attribute< Group, Element, TVR, VM::VM3_3n >	190
gdcm::Attribute< Group, Element, TVR, VM::VM3_n >	191
gdcm::AudioCodec	
AudioCodec	193
gdcm::Base64	
Class for Base64	195
gdcm::network::BaseCompositeMessage	
BaseCompositeMessage The Composite events described in section 3.7-2009 of the DICOM standard all use their own messages. These messages are constructed using Presentation Data Values, from section 3.8-2009 of the standard, and then fill in appropriate values in their datasets	196
gdcm::network::BasePDU	
BasePDU base class for PDUs	198
gdcm::BaseRootQuery	
BaseRootQuery contains: a baseclass which will produce a dataset for c-find and c-move with patient/study root	200
gdcm::SegmentHelper::BasicCodedEntry	
This structure defines a basic coded entry with all of its attributes	204
gdcm::BasicOffsetTable	
Class to represent a BasicOffsetTable	207
gdcm::Bitmap	
Bitmap class A bitmap based image. Used as parent for both IconImage and the main Pixel Data Image It does not contains any World Space information (IPP, IOP)	209
gdcm::BitmapToBitmapFilter	
BitmapToBitmapFilter class Super class for all filter taking an image and producing an output image	217
gdcm::BoxRegion	
Class for manipulation box region This is a very simple implementation of the Region class. It only support 3D box type region. It assumes the 3D Box does not have a tilt Origin is as (0,0,0)	219
gdcm::ByteBuffer	
ByteBuffer	222
gdcm::ByteSwap< T >	
ByteSwap	223
gdcm::ByteSwapFilter	
ByteSwapFilter In place byte-swapping of a dataset FIXME: FL status ??	224
gdcm::ByteValue	
Class to represent binary value (array of bytes)	224
gdcm::CAPICryptoFactory	229
gdcm::CAPICryptographicMessageSyntax	230
gdcm::network::CEchoRQ	
CEchoRQ this file defines the messages for the cecho action	233
gdcm::network::CEchoRSP	
CEchoRSP this file defines the messages for the cecho action	234
gdcm::network::CFind	235
gdcm::network::CFindCancelRQ	
CFindCancelRQ this file defines the messages for the cfind action	236
gdcm::network::CFindRQ	
CFindRQ this file defines the messages for the cfind action	237
gdcm::network::CFindRSP	
CFindRSP this file defines the messages for the cfind action	238

<a href="#">gdcm::network::CMoveCancelRq</a>	239
<a href="#">gdcm::network::CMoveRQ</a>	
CMoveRQ this file defines the messages for the cmove action	241
<a href="#">gdcm::network::CMoveRSP</a>	
CMoveRSP this file defines the messages for the cmove action	242
<a href="#">gdcm::Codec</a>	
Codec class	243
<a href="#">gdcm::Coder</a>	
Coder	244
<a href="#">gdcm::CodeString</a>	
CodeString This is an implementation of DICOM VR: CS The cstor will properly Trim so that operator== is correct	245
<a href="#">gdcm::Command</a>	
Command superclass for callback/observer methods	248
<a href="#">gdcm::CommandDataSet</a>	
Class to represent a <a href="#">Command DataSet</a>	250
<a href="#">gdcm::network::CompositeMessageFactory</a>	
CompositeMessageFactory This class constructs PDataPDUs, but that have been specifically constructed for the composite DICOM services (C-Echo, C-Find, C-Get, C-Move, and C-Store). It will also handle parsing the incoming data to determine which of the CompositePDUs the incoming data is, and so therefore allowing the scu to determine what to do with incoming data (if acting as a storescp server, for instance)	252
<a href="#">gdcm::CompositeNetworkFunctions</a>	
Composite Network Functions These functions provide a generic API to the DICOM functions implemented in GDCM. Advanced users can use this code as a template for building their own versions of these functions (for instance, to provide progress bars or some other way of handling returned query information), but for most users, these functions should be sufficient to interface with a PACS to a local machine. Note that these functions are not contained within a static class or some other class-style interface, because multiple connections can be instantiated in the same program. The DICOM standard is much more function oriented rather than class oriented in this instance, so the design of this API reflects that functional approach. These functions implements the following SCU operations:	253
<a href="#">gdcm::ConstCharWrapper</a>	
Do not use me	258
<a href="#">gdcm::CP246ExplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as CP246Explicit Data <a href="#">Element</a>	258
<a href="#">gdcm::CryptoFactory</a>	
Class to do handle the crypto factory	260
<a href="#">gdcm::CryptographicMessageSyntax</a>	262
<a href="#">gdcm::CSAElement</a>	
Class to represent a CSA <a href="#">Element</a>	264
<a href="#">gdcm::CSAHeader</a>	
Class for <a href="#">CSAHeader</a>	268
<a href="#">gdcm::CSAHeaderDict</a>	
Class to represent a map of <a href="#">CSAHeaderDictEntry</a>	272
<a href="#">gdcm::CSAHeaderDictEntry</a>	
Class to represent an Entry in the <a href="#">Dict</a> Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from <a href="#">gdcm::Tag</a> to the needed information	274
<a href="#">gdcm::CSAHeaderDictException</a>	276
<a href="#">gdcm::network::CStoreRQ</a>	
CStoreRQ this file defines the messages for the cecho action	277
<a href="#">gdcm::network::CStoreRSP</a>	
CStoreRSP this file defines the messages for the cecho action	278

<a href="#">gdcmm::Curve</a>	
<a href="#">Curve</a> class to handle element 50xx,3000 <a href="#">Curve</a> Data WARNING: This is deprecated and lastly defined in PS 3.3 - 2004	280
<a href="#">gdcmm::DataElement</a>	
Class to represent a Data <a href="#">Element</a> either Implicit or Explicit	282
<a href="#">gdcmm::DataElementException</a>	291
<a href="#">gdcmm::DataEvent</a>	
<a href="#">DataEvent</a>	292
<a href="#">gdcmm::DataSet</a>	
Class to represent a Data Set (which contains Data Elements) A Data Set represents an instance of a real world Information <a href="#">Object</a>	294
<a href="#">gdcmm::DataSetEvent</a>	
<a href="#">DataSetEvent</a> Special type of event triggered during the <a href="#">DataSet</a> store/move process	301
<a href="#">gdcmm::DataSetHelper</a>	
<a href="#">DataSetHelper</a> (internal class, not intended for user level)	304
<a href="#">gdcmm::Decoder</a>	
<a href="#">Decoder</a>	304
<a href="#">gdcmm::DefinedTerms</a>	
Defined Terms are used when the specified explicit Values may be extended by implementors to include additional new Values. These new Values shall be specified in the Conformance Statement (see PS 3.2) and shall not have the same meaning as currently defined Values in this standard. A Data <a href="#">Element</a> with Defined Terms that does not contain a <a href="#">Value</a> equivalent to one of the Values currently specified in this standard shall not be considered to have an invalid value. Note: Interpretation <a href="#">Type</a> ID (4008,0210) is an example of a Data <a href="#">Element</a> having Defined Terms. It is defined to have a <a href="#">Value</a> that may be one of the set of standard Values; REPORT or AMENDMENT (see PS 3.3). Because this Data <a href="#">Element</a> has Defined Terms other Interpretation <a href="#">Type</a> IDs may be defined by the implementor	306
<a href="#">gdcmm::Defs</a>	
FIXME I do not like the name ' <a href="#">Defs</a> '	306
<a href="#">gdcmm::DeltaEncodingCodec</a>	
<a href="#">DeltaEncodingCodec</a> compression used by some private vendor	309
<a href="#">gdcmm::DICOMDIR</a>	
<a href="#">DICOMDIR</a> class	310
<a href="#">gdcmm::DICOMDIRGenerator</a>	
<a href="#">DICOMDIRGenerator</a> class This is a STD-GEN-CD <a href="#">DICOMDIR</a> generator. ref: PS 3.11-2008 Annex D (Normative) - General Purpose CD-R and DVD Interchange Profiles	311
<a href="#">gdcmm::Dict</a>	
Class to represent a map of <a href="#">DictEntry</a>	313
<a href="#">gdcmm::DictConverter</a>	
Class to convert a .dic file into something else:	316
<a href="#">gdcmm::DictEntry</a>	
Class to represent an Entry in the <a href="#">Dict</a> Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from <a href="#">gdcmm::Tag</a> to the needed information	318
<a href="#">gdcmm::DictPrinter</a>	
<a href="#">DictPrinter</a> class	320
<a href="#">gdcmm::Dicts</a>	
Class to manipulate the sum of knowledge (all the dict user load)	322
<a href="#">gdcmm::network::DIMSE</a>	
<a href="#">DIMSE</a> PS 3.7 - 2009 Annex E <a href="#">Command</a> Dictionary (Normative) E.1 REGISTRY OF DICOM CO↔MMAND ELEMENTS <a href="#">Table</a> E.1-1 COMMAND FIELDS (PART 1)	325
<a href="#">gdcmm::DirectionCosines</a>	
Class to handle <a href="#">DirectionCosines</a>	326
<a href="#">gdcmm::Directory</a>	
Class for manipulation directories	328

<a href="#">gdcm::DirectoryHelper</a>	
<a href="#">DirectoryHelper</a>	this class is designed to help mitigate some of the commonly performed operations on directories. namely: 1) the ability to determine the number of series in a directory by what type of series is present 2) the ability to find all ct series in a directory 3) the ability to find all mr series in a directory 4) to load a set of DataSets from a series that's already been sorted by the IPP sorter 5) For rtstruct stuff, you need to know the sopinstanceuid of each z plane, so there's a retrieval function for that 6) then a few other functions for rtstruct writeouts . . . . .
	331
<a href="#">gdcm::DummyValueGenerator</a>	
Class for generating dummy value . . . . .	332
<a href="#">gdcm::Dumper</a>	
Codec class . . . . .	333
<a href="#">gdcm::Element&lt; TVR, TVM &gt;</a>	
Element class . . . . .	335
<a href="#">gdcm::Element&lt; TVR, VM::VM1_2 &gt;</a>	338
<a href="#">gdcm::Element&lt; TVR, VM::VM1_n &gt;</a>	339
<a href="#">gdcm::Element&lt; TVR, VM::VM2_2n &gt;</a>	341
<a href="#">gdcm::Element&lt; TVR, VM::VM2_n &gt;</a>	343
<a href="#">gdcm::Element&lt; TVR, VM::VM3_3n &gt;</a>	344
<a href="#">gdcm::Element&lt; TVR, VM::VM3_n &gt;</a>	346
<a href="#">gdcm::Element&lt; VR::AS, VM::VM5 &gt;</a>	347
<a href="#">gdcm::Element&lt; VR::OB, VM::VM1 &gt;</a>	348
<a href="#">gdcm::Element&lt; VR::OW, VM::VM1 &gt;</a>	349
<a href="#">gdcm::ElementDisableCombinations&lt; TVR, TVM &gt;</a>	
A class which is used to produce compile errors for an invalid combination of template parameters . . . . .	351
<a href="#">gdcm::ElementDisableCombinations&lt; VR::OB, VM::VM1_n &gt;</a>	352
<a href="#">gdcm::ElementDisableCombinations&lt; VR::OW, VM::VM1_n &gt;</a>	352
<a href="#">gdcm::EncapsulatedDocument</a>	
EncapsulatedDocument . . . . .	352
<a href="#">gdcm::EncodingImplementation&lt; T &gt;</a>	
EncodingImplementation . . . . .	353
<a href="#">gdcm::EncodingImplementation&lt; VR::VRASCII &gt;</a>	353
<a href="#">gdcm::EncodingImplementation&lt; VR::VRBINARY &gt;</a>	354
<a href="#">gdcm::EndEvent</a>	355
<a href="#">gdcm::EnumeratedValues</a>	
Element. A Data <a href="#">Element</a> with Enumerated Values that does not have a <a href="#">Value</a> equivalent to one of the Values specified in this standard has an invalid value within the scope of a specific Information Object/SOP Class definition. Note: . . . . .	356
<a href="#">gdcm::Event</a>	
Superclass for callback/observer methods . . . . .	357
<a href="#">gdcm::Exception</a>	
Exception . . . . .	359
<a href="#">gdcm::ExitEvent</a>	361
<a href="#">gdcm::ExplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as Explicit Data <a href="#">Element</a> . . . . .	362
<a href="#">gdcm::ExplicitImplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as ExplicitImplicit Data <a href="#">Element</a> . . . . .	363
<a href="#">gdcm::Fiducials</a>	
Fiducials . . . . .	365
<a href="#">gdcm::File</a>	
DICOM <a href="#">File</a> See PS 3.10 <a href="#">File</a> : A <a href="#">File</a> is an ordered string of zero or more bytes, where the first byte is at the beginning of the file and the last byte at the end of the <a href="#">File</a> . Files are identified by a unique <a href="#">File</a> ID and may be written, read and/or deleted . . . . .	366
<a href="#">gdcm::FileAnonymizer</a>	
FileAnonymizer . . . . .	369

<a href="#">gdcm::FileChangeTransferSyntax</a>	
<a href="#">FileChangeTransferSyntax</a> . . . . .	372
<a href="#">gdcm::FileDerivation</a>	
<a href="#">FileDerivation</a> class See PS 3.16 - 2008 For the list of Code <a href="#">Value</a> that can be used for in Derivation	
Code Sequence . . . . .	374
<a href="#">gdcm::FileExplicitFilter</a>	
<a href="#">FileExplicitFilter</a> class After changing a file from Implicit to Explicit representation (see <a href="#">Image↔ChangeTransferSyntax</a> ) one operation is to make sure the <a href="#">VR</a> of each DICOM attribute are accurate and do match the one from PS 3.6. Indeed when a file is written in Implicit representation, the <a href="#">VR</a> is not stored directly in the file . . . . .	377
<a href="#">gdcm::FileMetaInformation</a>	
Class to represent a <a href="#">File</a> Meta Information . . . . .	379
<a href="#">gdcm::Filename</a>	
Class to manipulate file name's . . . . .	384
<a href="#">gdcm::FileNameEvent</a>	
<a href="#">FileNameEvent</a> Special type of event triggered during processing of <a href="#">FileSet</a> . . . . .	386
<a href="#">gdcm::FilenameGenerator</a>	
<a href="#">FilenameGenerator</a> . . . . .	388
<a href="#">gdcm::FileSet</a>	
File-set: A File-set is a collection of DICOM Files (and possibly non-DICOM Files) that share a common naming space within which <a href="#">File</a> IDs are unique . . . . .	391
<a href="#">gdcm::FileStreamer</a>	
<a href="#">FileStreamer</a> This class let a user create a massive DICOM <a href="#">DataSet</a> from a template DICOM file, by appending chunks of data . . . . .	392
<a href="#">gdcm::FileWithName</a>	
<a href="#">FileWithName</a> . . . . .	396
<a href="#">gdcm::FindPatientRootQuery</a>	
<a href="#">PatientRootQuery</a> contains: the class which will produce a dataset for c-find with patient root . . . .	397
<a href="#">gdcm::FindStudyRootQuery</a>	
<a href="#">FindStudyRootQuery</a> contains: the class which will produce a dataset for C-FIND with study root . .	399
<a href="#">gdcm::Fragment</a>	
Class to represent a <a href="#">Fragment</a> . . . . .	401
<a href="#">gdcm::Global</a>	
<a href="#">Global</a> . . . . .	404
<a href="#">gdcm::GroupDict</a>	
Class to represent the mapping from group number to its abbreviation and name . . . . .	406
<a href="#">gdcm::IconImageFilter</a>	
<a href="#">IconImageFilter</a> This filter will extract icons from a <a href="#">File</a> This filter will loop over all known sequence (public and private) that may contains an <a href="#">IconImage</a> and retrieve them. The filter will fails with a value of false if no icon can be found Since it handle both public and private icon type, one should not assume the icon is in uncompress form, some private vendor store private icon in JPEG8/JPEG12	408
<a href="#">gdcm::IconImageGenerator</a>	
<a href="#">IconImageGenerator</a> This filter will generate a valid Icon from the Pixel Data element (an instance of <a href="#">Pixmap</a> ). To generate a valid Icon, one is only allowed the following Photometric Interpretation: . .	410
<a href="#">gdcm::ignore_char</a> . . . . .	412
<a href="#">gdcm::Image</a>	
<a href="#">Image</a> This is the container for an <a href="#">Image</a> in the general sense. From this container you should be able to request information like: . . . . .	413
<a href="#">gdcm::ImageApplyLookupTable</a>	
<a href="#">ImageApplyLookupTable</a> class It applies the LUT the <a href="#">PixelData</a> (only PALETTE_COLOR images) Output will be a <a href="#">PhotometricInterpretation</a> =RGB image . . . . .	417
<a href="#">gdcm::ImageChangePhotometricInterpretation</a>	
<a href="#">ImageChangePhotometricInterpretation</a> class Class to change the Photometric Interpretation of an input DICOM . . . . .	419



gdcm::ImageChangePlanarConfiguration	
ImageChangePlanarConfiguration	class Class to change the Planar configuration of an input DICOM
By default it will change into the more usual representation: PlanarConfiguration = 0	422
gdcm::ImageChangeTransferSyntax	
ImageChangeTransferSyntax	class Class to change the transfer syntax of an input DICOM
	425
gdcm::ImageCodec	
ImageCodec	429
gdcm::ImageConverter	
Image Converter	435
gdcm::ImageFragmentSplitter	
ImageFragmentSplitter	class For single frame image, DICOM standard allow splitting the frame into multiple fragments
	436
gdcm::ImageHelper	
ImageHelper	(internal class, not intended for user level)
	439
gdcm::ImageReader	
ImageReader	443
gdcm::ImageRegionReader	
ImageRegionReader	446
gdcm::ImageToImageFilter	
ImageToImageFilter	class Super class for all filter taking an image and producing an output image
	449
gdcm::ImageWriter	
ImageWriter	451
gdcm::network::ImplementationClassUIDSub	
ImplementationClassUIDSub	PS 3.7 Table D.3-1 IMPLEMENTATION CLASS UID SUB-ITEM FIELDS (A-ASSOCIATE-RQ)
	454
gdcm::network::ImplementationUIDSub	
ImplementationUIDSub	Table D.3-2 IMPLEMENTATION UID SUB-ITEM FIELDS (A-ASSOCIATE-RQ)
	454
gdcm::network::ImplementationVersionNameSub	
ImplementationVersionNameSub	Table D.3-3 IMPLEMENTATION VERSION NAME SUB-ITEM FIELDS (A-ASSOCIATE-RQ)
	455
gdcm::ImplicitDataElement	
Class to represent an <i>Implicit VR</i> Data Element	456
gdcm::InitializeEvent	458
gdcm::IOD	
Class for representing a IOD	459
gdcm::IODEntry	
Class for representing a IODEntry	460
gdcm::IODs	
Class for representing a IODs	462
gdcm::IPPSorter	
IPPSorter	Implement a simple Image Position (Patient) sorter, along the Image Orientation (Patient) direction. This algorithm does NOT support duplicate and will FAIL in case of duplicate IPP
	464
gdcm::Item	
Class to represent an Item	A component of the value of a Data Element that is of Value Representation Sequence of Items. An Item contains a Data Set. See PS 3.5 7.5.1 Item Encoding Rules
Each Item of a Data Element of VR SQ shall be encoded as a DICOM Standard Data Element with a specific Data Element Tag of Value (FFFE,E000). The Item Tag is followed by a 4 byte Item Length field encoded in one of the following two ways Explicit/ Implicit	467
gdcm::IterationEvent	470
gdcm::JPEG12Codec	
Class to do JPEG 12bits (lossy & lossless)	471
gdcm::JPEG16Codec	
Class to do JPEG 16bits (lossless)	474

<a href="#">gdcm::JPEG2000Codec</a>	
Class to do JPEG 2000	476
<a href="#">gdcm::JPEG8Codec</a>	
Class to do JPEG 8bits (lossy & lossless)	479
<a href="#">gdcm::JPEGCodec</a>	
JPEG codec Class to do JPEG (8bits, 12bits, 16bits lossy & lossless). It redispach in between the different codec implementation: <a href="#">JPEG8Codec</a> , <a href="#">JPEG12Codec</a> & <a href="#">JPEG16Codec</a> It also support inconsistency in between DICOM header and JPEG compressed stream <a href="#">ImageCodec</a> implementation for the JPEG case	482
<a href="#">gdcm::JPEGLSCodec</a>	
JPEG-LS	486
<a href="#">gdcm::JSON</a>	490
<a href="#">gdcm::KAKADUCodec</a>	
<a href="#">KAKADUCodec</a>	491
<a href="#">gdcm::LO</a>	
LO	493
<a href="#">gdcm::LookupTable</a>	
<a href="#">LookupTable</a> class	496
<a href="#">gdcm::Scanner::ltstr</a>	500
<a href="#">gdcm::Macro</a>	
Class for representing a <a href="#">Macro</a>	501
<a href="#">gdcm::Macros</a>	
Class for representing a <a href="#">Modules</a>	502
<a href="#">gdcm::network::MaximumLengthSub</a>	
<a href="#">MaximumLengthSub</a> Annex D Table D.1-1 MAXIMUM LENGTH SUB-ITEM FIELDS (A-ASSOCIATION ← TE-RQ)	503
<a href="#">gdcm::MD5</a>	
Class for <a href="#">MD5</a>	504
<a href="#">gdcm::MediaStorage</a>	
<a href="#">MediaStorage</a>	505
<a href="#">gdcm::MemberCommand&lt; T &gt;</a>	
<a href="#">Command</a> subclass that calls a pointer to a member function	512
<a href="#">gdcm::MeshPrimitive</a>	
This class defines surface mesh primitives. It is designed from surface mesh primitives macro	516
<a href="#">gdcm::ModifiedEvent</a>	519
<a href="#">gdcm::Module</a>	
Class for representing a <a href="#">Module</a>	521
<a href="#">gdcm::ModuleEntry</a>	
Class for representing a <a href="#">ModuleEntry</a>	522
<a href="#">gdcm::Modules</a>	
Class for representing a <a href="#">Modules</a>	525
<a href="#">gdcm::MovePatientRootQuery</a>	
<a href="#">MovePatientRootQuery</a> contains: the class which will produce a dataset for c-move with patient root	527
<a href="#">gdcm::MoveStudyRootQuery</a>	
<a href="#">MoveStudyRootQuery</a> contains: the class which will produce a dataset for C-MOVE with study root	529
<a href="#">gdcm::NestedModuleEntries</a>	
Class for representing a <a href="#">NestedModuleEntries</a>	531
<a href="#">gdcm::NoEvent</a>	533
<a href="#">gdcm::Object</a>	
<a href="#">Object</a>	534
<a href="#">gdcm::OpenSSLCryptoFactory</a>	537
<a href="#">gdcm::OpenSSLCryptographicMessageSyntax</a>	538
<a href="#">gdcm::OpenSSLP7CryptoFactory</a>	540

<a href="#">gdcm::OpenSSLP7CryptographicMessageSyntax</a>	
Class for <a href="#">CryptographicMessageSyntax</a> encryption. This is just a simple wrapper around openssl PKCS7_encrypt functionalities	542
<a href="#">gdcm::Orientation</a>	
Class to handle <a href="#">Orientation</a>	544
<a href="#">gdcm::Overlay</a>	
Overlay class	546
<a href="#">gdcm::ParseException</a>	
<a href="#">ParseException</a> Standard exception handling object	552
<a href="#">gdcm::Parser</a>	
Parser ala XML_Parser from expat (SAX)	554
<a href="#">gdcm::Patient</a>	
See PS 3.3 - 2007 DICOM MODEL OF THE REAL-WORLD, p 54	556
<a href="#">gdcm::network::PDataTFPDU</a>	
<a href="#">PDataTFPDU</a> Table 9-22 P-DATA-TF PDU FIELDS	556
<a href="#">gdcm::PDBElement</a>	
Class to represent a PDB <a href="#">Element</a>	559
<a href="#">gdcm::PDBHeader</a>	
Class for <a href="#">PDBHeader</a>	561
<a href="#">gdcm::PDFCodec</a>	
<a href="#">PDFCodec</a> class	562
<a href="#">gdcm::network::PDUFactory</a>	
<a href="#">PDUFactory</a> basically, given an initial byte, construct the appropriate PDU. This way, the event loop doesn't have to know about all the different PDU types	564
<a href="#">gdcm::PersonName</a>	
<a href="#">PersonName</a> class	566
<a href="#">gdcm::PGXCodec</a>	
Class to do PGX See PGX as used in JPEG 2000 implementation and reference images	567
<a href="#">gdcm::PhotometricInterpretation</a>	
Class to represent an <a href="#">PhotometricInterpretation</a>	569
<a href="#">gdcm::PixelFormat</a>	
<a href="#">PixelFormat</a>	571
<a href="#">gdcm::Pixmap</a>	
<a href="#">Pixmap</a> class A bitmap based image. Used as parent for both <a href="#">IconImage</a> and the main <a href="#">Pixel Data Image</a> It does not contains any World Space information (IPP, IOP)	577
<a href="#">gdcm::PixmapReader</a>	
<a href="#">PixmapReader</a>	579
<a href="#">gdcm::PixmapToPixmapFilter</a>	
<a href="#">PixmapToPixmapFilter</a> class Super class for all filter taking an image and producing an output image	583
<a href="#">gdcm::PixmapWriter</a>	
<a href="#">PixmapWriter</a> This class will takes two inputs:	585
<a href="#">gdcm::PNMCodec</a>	
Class to do PNM PNM is the Portable anymap file format. The main web page can be found at: <a href="http://netpbm.sourceforge.net/">http://netpbm.sourceforge.net/</a>	588
<a href="#">gdcm::Preamble</a>	
DICOM <a href="#">Preamble</a> (Part 10)	591
<a href="#">gdcm::PresentationContext</a>	
<a href="#">PresentationContext</a>	592
<a href="#">gdcm::network::PresentationContextAC</a>	
<a href="#">PresentationContextAC</a> Table 9-18 PRESENTATION CONTEXT ITEM FIELDS	594

<a href="#">gdcm::PresentationContextGenerator</a>	
<a href="#">PresentationContextGenerator</a>	This class is responsible for generating the proper <a href="#">Presentation↔Context</a> that will be used in subsequent operation during a DICOM Query/Retrieve association. The step of the association is very sensible as special care need to be taken to explicitly define what instance are going to be send and how they are encoded . . . . .
	595
<a href="#">gdcm::network::PresentationContextRQ</a>	
<a href="#">PresentationContextRQ</a>	<a href="#">Table 9-13 PRESENTATION CONTEXT ITEM FIELDS</a> . . . . .
	597
<a href="#">gdcm::network::PresentationDataValue</a>	
<a href="#">PresentationDataValue</a>	<a href="#">Table 9-23 PRESENTATION-DATA-VALUE ITEM FIELDS</a> . . . . .
	599
<a href="#">gdcm::Printer</a>	
<a href="#">Printer</a>	class . . . . .
	601
<a href="#">gdcm::PrivateDict</a>	
<a href="#">Private</a>	<a href="#">Dict</a> . . . . .
	604
<a href="#">gdcm::PrivateTag</a>	
	Class to represent a Private DICOM Data <a href="#">Element</a> ( <a href="#">Attribute</a> ) <a href="#">Tag</a> (Group, <a href="#">Element</a> , Owner) . . . . .
	606
<a href="#">gdcm::ProgressEvent</a>	
<a href="#">ProgressEvent</a>	Special type of event triggered during . . . . .
	608
<a href="#">gdcm::PVRGCodec</a>	
<a href="#">PVRGCodec</a>	. . . . .
	610
<a href="#">gdcm::PythonFilter</a>	
<a href="#">PythonFilter</a>	<a href="#">PythonFilter</a> is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a <a href="#">DataElement</a> into a string, typically this is a nice feature to have for wrapped language . . . . .
	612
<a href="#">gdcm::QueryBase</a>	
<a href="#">QueryBase</a>	contains: the base class for constructing a query dataset for a C-FIND and a C-MOVE . . . . .
	613
<a href="#">gdcm::QueryFactory</a>	
<a href="#">QueryFactory.h</a>	. . . . .
	615
<a href="#">gdcm::QueryImage</a>	
<a href="#">QueryImage</a>	contains: class to construct an image-based query for C-FIND and C-MOVE . . . . .
	616
<a href="#">gdcm::QueryPatient</a>	
<a href="#">QueryPatient</a>	contains: class to construct a patient-based query for c-find and c-move . . . . .
	618
<a href="#">gdcm::QuerySeries</a>	
<a href="#">QuerySeries</a>	contains: class to construct a series-based query for c-find and c-move . . . . .
	620
<a href="#">gdcm::QueryStudy</a>	
<a href="#">QueryStudy.h</a>	contains: class to construct a study-based query for C-FIND and C-MOVE . . . . .
	622
<a href="#">gdcm::RAWCodec</a>	
<a href="#">RAWCodec</a>	class . . . . .
	624
<a href="#">gdcm::Reader</a>	
<a href="#">Reader</a>	ala DOM (Document <a href="#">Object</a> Model) . . . . .
	627
<a href="#">gdcm::Region</a>	
	Class for manipulation region . . . . .
	632
<a href="#">gdcm::Rescaler</a>	
<a href="#">Rescale</a>	class This class is meant to apply the linear transform of Stored Pixel <a href="#">Value</a> to Real World <a href="#">Value</a> . This is mostly found in CT or PET dataset, where the value are stored using one type, but need to be converted to another scale using a linear transform. There are basically two cases: In CT: the linear transform is generally integer based. E.g. the Stored Pixel <a href="#">Type</a> is unsigned short 12bits, but to get Hounsfield unit, one need to apply the linear transform:
	$RWV = 1. * SV - 1024$
	So the best scalar to store the Real World <a href="#">Value</a> will be 16 bits signed type . . . . .
	634
<a href="#">gdcm::RLECodec</a>	
	Class to do RLE . . . . .
	637

<a href="#">gdcm::network::RoleSelectionSub</a>	
<a href="#">RoleSelectionSub</a> PS 3.7 Table D.3-9 SCP/SCU ROLE SELECTION SUB-ITEM FIELDS (A-ASSOCIATE-RQ)	640
<a href="#">gdcm::SerieHelper::Rule</a>	641
<a href="#">gdcm::Scanner</a>	
<a href="#">Scanner</a> This filter is meant for quickly browsing a <a href="#">FileSet</a> (a set of files on disk). Special consideration are taken so as to read the minimum amount of information in each file in order to retrieve the user specified set of DICOM <a href="#">Attribute</a>	642
<a href="#">gdcm::Segment</a>	
This class defines a segment. It mainly contains attributes of group 0x0062. In addition, it can be associated with surface	648
<a href="#">gdcm::SegmentedPaletteColorLookupTable</a>	
<a href="#">SegmentedPaletteColorLookupTable</a> class	652
<a href="#">gdcm::SegmentReader</a>	
This class defines a segment reader. It reads attributes of group 0x0062	654
<a href="#">gdcm::SegmentWriter</a>	
This class defines a segment writer. It writes attributes of group 0x0062	656
<a href="#">gdcm::SequenceOfFragments</a>	
Class to represent a Sequence Of Fragments	659
<a href="#">gdcm::SequenceOfItems</a>	
Class to represent a Sequence Of Items (value representation : SQ)	664
<a href="#">gdcm::SerieHelper</a>	
<a href="#">SerieHelper</a> DO NOT USE this class, it is only a temporary solution for ITK migration from GDCM 1.x to GDCM 2.x It will disappear soon, you've been warned	670
<a href="#">gdcm::Series</a>	
<a href="#">Series</a>	672
<a href="#">gdcm::network::ServiceClassApplicationInformation</a>	673
<a href="#">gdcm::ServiceClassUser</a>	
<a href="#">ServiceClassUser</a>	674
<a href="#">gdcm::SHA1</a>	
Class for <a href="#">SHA1</a>	679
<a href="#">gdcm::SimpleMemberCommand&lt; T &gt;</a>	
<a href="#">Command</a> subclass that calls a pointer to a member function	680
<a href="#">gdcm::SimpleSubjectWatcher</a>	
<a href="#">SimpleSubjectWatcher</a> This is a typical <a href="#">Subject</a> Watcher class. It will observe all events	684
<a href="#">gdcm::SmartPointer&lt; ObjectType &gt;</a>	
Class for Smart Pointer	685
<a href="#">gdcm::network::SOPClassExtendedNegociationSub</a>	
<a href="#">SOPClassExtendedNegociationSub</a> PS 3.7 Table D.3-11 SOP CLASS EXTENDED NEGOTIATION SUB-ITEM FIELDS (A-ASSOCIATE-RQ and A-ASSOCIATE-AC)	688
<a href="#">gdcm::SOPClassUIDToIOD</a>	
Class convert a class SOP Class UID into <a href="#">IOD</a>	689
<a href="#">gdcm::Sorter</a>	
<a href="#">Sorter</a> General class to do sorting using a custom function You simply need to provide a function of type: <a href="#">Sorter::SortFunction</a>	690
<a href="#">gdcm::Spacing</a>	
Class for <a href="#">Spacing</a>	694
<a href="#">gdcm::Spectroscopy</a>	
<a href="#">Spectroscopy</a> class	696
<a href="#">gdcm::SplitMosaicFilter</a>	
<a href="#">SplitMosaicFilter</a> class Class to reshuffle bytes for a SIEMENS Mosaic image Siemens CSA <a href="#">Image</a> Header CSA:= Common Siemens Architecture, sometimes also known as Common syngo Architecture	696
<a href="#">gdcm::StartEvent</a>	697

<a href="#">gdcmm::static_assert_test&lt; x &gt;</a>	699
<a href="#">gdcmm::STATIC_ASSERTION_FAILURE&lt; x &gt;</a>	699
<a href="#">gdcmm::STATIC_ASSERTION_FAILURE&lt; true &gt;</a>	699
<a href="#">gdcmm::StreamImageReader</a>	
<a href="#">StreamImageReader</a>	699
<a href="#">gdcmm::StreamImageWriter</a>	
<a href="#">StreamImageReader</a>	702
<a href="#">gdcmm::String&lt; TDelimiter, TMaxLength, TPadChar &gt;</a>	
<a href="#">String</a>	707
<a href="#">gdcmm::StringFilter</a>	
<a href="#">StringFilter</a> <a href="#">StringFilter</a> is the class that make gdcmm2.x looks more like gdcmm1 and transform the binary blob contained in a <a href="#">DataElement</a> into a string, typically this is a nice feature to have for wrapped language	710
<a href="#">gdcmm::Study</a>	
<a href="#">Study</a>	713
<a href="#">gdcmm::Subject</a>	
<a href="#">Subject</a>	713
<a href="#">gdcmm::Surface</a>	
This class defines a SURFACE IE. This members are taken from required surface mesh module attributes	716
<a href="#">gdcmm::SurfaceHelper</a>	
<a href="#">SurfaceHelper</a> Helper class for <a href="#">Surface</a> object	722
<a href="#">gdcmm::SurfaceReader</a>	
This class defines a SURFACE IE reader. It reads surface mesh module attributes	725
<a href="#">gdcmm::SurfaceWriter</a>	
This class defines a SURFACE IE writer. It writes surface mesh module attributes	727
<a href="#">gdcmm::SwapCode</a>	
<a href="#">SwapCode</a> representation	729
<a href="#">gdcmm::SwapperDoOp</a>	731
<a href="#">gdcmm::SwapperNoOp</a>	732
<a href="#">gdcmm::System</a>	
Class to do system operation	732
<a href="#">gdcmm::Table</a>	
<a href="#">Table</a>	736
<a href="#">gdcmm::TableEntry</a>	
<a href="#">TableEntry</a>	737
<a href="#">gdcmm::TableReader</a>	
Class for representing a <a href="#">TableReader</a>	738
<a href="#">gdcmm::network::TableRow</a>	740
<a href="#">gdcmm::Tag</a>	
Class to represent a DICOM Data <a href="#">Element</a> ( <a href="#">Attribute</a> ) <a href="#">Tag</a> (Group, <a href="#">Element</a> ). Basically an uint32_t which can also be expressed as two uint16_t (group and element)	741
<a href="#">gdcmm::TagPath</a>	
Class to handle a path of tag	748
<a href="#">gdcmm::Testing</a>	
Class for testing	749
<a href="#">gdcmm::Trace</a>	
<a href="#">Trace</a>	754
<a href="#">gdcmm::TransferSyntax</a>	
Class to manipulate Transfer Syntax	757
<a href="#">gdcmm::network::TransferSyntaxSub</a>	
<a href="#">TransferSyntaxSub</a> <a href="#">Table</a> 9-15 TRANSFER SYNTAX SUB-ITEM FIELDS	761
<a href="#">gdcmm::network::Transition</a>	762

gdcmm::Type	
Type	764
gdcmm::UI	765
gdcmm::UIDGenerator	
Class for generating unique UID	766
gdcmm::UIDs	
All known uids	768
gdcmm::network::ULAction	
ULAction A <a href="#">ULConnection</a> in a given ULState can perform certain ULActions. This base class provides the interface for running those ULActions on a given <a href="#">ULConnection</a>	786
gdcmm::network::ULActionAA1	789
gdcmm::network::ULActionAA2	790
gdcmm::network::ULActionAA3	791
gdcmm::network::ULActionAA4	792
gdcmm::network::ULActionAA5	793
gdcmm::network::ULActionAA6	794
gdcmm::network::ULActionAA7	795
gdcmm::network::ULActionAA8	796
gdcmm::network::ULActionAE1	797
gdcmm::network::ULActionAE2	798
gdcmm::network::ULActionAE3	799
gdcmm::network::ULActionAE4	800
gdcmm::network::ULActionAE5	801
gdcmm::network::ULActionAE6	802
gdcmm::network::ULActionAE7	803
gdcmm::network::ULActionAE8	804
gdcmm::network::ULActionAR1	805
gdcmm::network::ULActionAR10	806
gdcmm::network::ULActionAR2	807
gdcmm::network::ULActionAR3	808
gdcmm::network::ULActionAR4	809
gdcmm::network::ULActionAR5	810
gdcmm::network::ULActionAR6	811
gdcmm::network::ULActionAR7	812
gdcmm::network::ULActionAR8	813
gdcmm::network::ULActionAR9	814
gdcmm::network::ULActionDT1	815
gdcmm::network::ULActionDT2	816
gdcmm::network::ULBasicCallback	
ULBasicCallback This is the most basic of callbacks for how the <a href="#">ULConnectionManager</a> handles incoming datasets. DataSets are just concatenated to the mDataSets vector, and the result can be pulled out of the vector by later code. Alternatives to this method include progress updates, saving to disk, etc. This class is NOT THREAD SAFE. Access the dataset vector after the entire set of datasets has been returned by the <a href="#">ULConnectionManager</a>	817
gdcmm::network::ULConnection	
ULConnection This is the class that contains the socket to another machine, and passes data through itself, as well as maintaining a sense of state	818
gdcmm::network::ULConnectionCallback	821
gdcmm::network::ULConnectionInfo	
ULConnectionInfo this class contains all the information about a particular connection as established by the user. That is, it's: User Information Calling AE Title Called AE Title IP address/computer name IP Port A connection must be established with this information, that's subsequently placed into various primitives for actual communication	822

<a href="#">gdcm::network::ULConnectionManager</a>	
<a href="#">ULConnectionManager</a>	The <a href="#">ULConnectionManager</a> performs actions on the <a href="#">ULConnection</a> given inputs from the user and from the state of what's going on around the connection (ie, timeouts of the ARTIM timer, responses from the peer across the connection, etc) . . . . .
	824
<a href="#">gdcm::network::ULEvent</a>	
<a href="#">ULEvent</a>	base class for network events . . . . .
	826
<a href="#">gdcm::network::ULTransitionTable</a>	
<a href="#">ULTransitionTable</a>	The transition table of all the ULEvents, new ULActions, and ULStates . . . . .
	827
<a href="#">gdcm::network::ULWritingCallback</a>	. . . . .
	828
<a href="#">gdcm::UNExplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as UNExplicit Data <a href="#">Element</a>	. . . . .
	830
<a href="#">gdcm::UNExplicitImplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as ExplicitImplicit Data <a href="#">Element</a>	This class gather two known bugs: . . . . .
	831
<a href="#">gdcm::Unpacker12Bits</a>	
Pack/Unpack 12 bits pixel into 16bits . . . . .	833
<a href="#">gdcm::Usage</a>	
Usage . . . . .	834
<a href="#">gdcm::UserEvent</a>	. . . . .
	836
<a href="#">gdcm::network::UserInformation</a>	
<a href="#">UserInformation Table</a>	9-16 USER INFORMATION ITEM FIELDS . . . . .
	837
<a href="#">gdcm::UUIDGenerator</a>	
Class for generating unique UUID generate DCE 1.1 uid . . . . .	838
<a href="#">gdcm::Validate</a>	
Validate class . . . . .	839
<a href="#">gdcm::Value</a>	
Class to represent the value of a Data <a href="#">Element</a> . . . . .	840
<a href="#">gdcm::ValueIO&lt; TDE, TSwap, TType &gt;</a>	
Class to dispatch template calls . . . . .	842
<a href="#">gdcm::Version</a>	
Major/minor and build version . . . . .	843
<a href="#">gdcm::VL</a>	
Value Length . . . . .	844
<a href="#">gdcm::VM</a>	
Value Multiplicity Looking at the DICOMV3 dict only there is very few cases: 1 2 3 4 5 6 8 16 24 1-2 1-3 1-8 1-32 1-99 1-n 2-2n 2-n 3-3n 3-n . . . . .	846
<a href="#">gdcm::VMToLength&lt; T &gt;</a>	. . . . .
	850
<a href="#">gdcm::VR</a>	
VR class This is adapted from DICOM standard The biggest difference is the INVALID VR and the composite one that differ from standard (more like an addition) This allow us to represent all the possible case express in the DICOMV3 dict . . . . .	850
<a href="#">gdcm::VR16ExplicitDataElement</a>	
Class to read/write a <a href="#">DataElement</a> as Explicit Data <a href="#">Element</a> . . . . .	855
<a href="#">gdcm::VRToEncoding&lt; T &gt;</a>	. . . . .
	857
<a href="#">gdcm::VRToType&lt; T &gt;</a>	. . . . .
	857
<a href="#">gdcm::VRVLSize&lt; T &gt;</a>	. . . . .
	858
<a href="#">gdcm::VRVLSize&lt; 0 &gt;</a>	. . . . .
	858
<a href="#">gdcm::VRVLSize&lt; 1 &gt;</a>	. . . . .
	858
<a href="#">vtkGDCMImageReader</a>	. . . . .
	859
<a href="#">vtkGDCMImageReader2</a>	. . . . .
	865
<a href="#">vtkGDCMImageWriter</a>	. . . . .
	871
<a href="#">vtkGDCMMedicalImageProperties</a>	. . . . .
	875
<a href="#">vtkGDCMPolyDataReader</a>	. . . . .
	877
<a href="#">vtkGDCMPolyDataWriter</a>	. . . . .
	880



<a href="#">vtkGDCMTesting</a>	883
<a href="#">vtkGDCMThreadedImageReader</a>	886
<a href="#">vtkGDCMThreadedImageReader2</a>	888
<a href="#">vtkImageColorViewer</a>	892
<a href="#">vtkImageMapToColors16</a>	899
<a href="#">vtkImageMapToWindowLevelColors2</a>	902
<a href="#">vtkImagePlanarComponentsToComponents</a>	904
<a href="#">vtkImageRGBToYBR</a>	906
<a href="#">vtkImageYBRToRGB</a>	908
<a href="#">vtkLookupTable16</a>	909
<a href="#">vtkRTStructSetProperties</a>	911
<a href="#">gdcm::Waveform</a>	
<a href="#">Waveform</a> class	916
<a href="#">gdcm::Writer</a>	
<a href="#">Writer</a> ala DOM (Document <a href="#">Object</a> Model) This class is a non-validating writer, it will only performs well- formedness check only	916
<a href="#">gdcm::XMLDictReader</a>	
Class for representing a <a href="#">XMLDictReader</a>	921
<a href="#">gdcm::XMLPrinter</a>	923
<a href="#">gdcm::XMLPrivateDictReader</a>	
Class for representing a <a href="#">XMLPrivateDictReader</a>	926



## Chapter 25

# File Index

### 25.1 File List

Here is a list of all files with brief descriptions:

<a href="#">gdcmAAbortPDU.h</a>	929
<a href="#">gdcmAAssociateACPDU.h</a>	930
<a href="#">gdcmAAssociateRJPDU.h</a>	930
<a href="#">gdcmAAssociateRQPDU.h</a>	931
<a href="#">gdcmAbstractSyntax.h</a>	932
<a href="#">gdcmAnonymizeEvent.h</a>	933
<a href="#">gdcmAnonymizer.h</a>	934
<a href="#">gdcmApplicationContext.h</a>	935
<a href="#">gdcmApplicationEntity.h</a>	936
<a href="#">gdcmAReleaseRPPDU.h</a>	936
<a href="#">gdcmAReleaseRQPDU.h</a>	937
<a href="#">gdcmARTIMTimer.h</a>	938
<a href="#">gdcmASN1.h</a>	939
<a href="#">gdcmAsynchronousOperationsWindowSub.h</a>	939
<a href="#">gdcmAttribute.h</a>	940
<a href="#">gdcmAudioCodec.h</a>	942
<a href="#">gdcmBase64.h</a>	942
<a href="#">gdcmBaseCompositeMessage.h</a>	943
<a href="#">gdcmBasePDU.h</a>	944
<a href="#">gdcmBaseRootQuery.h</a>	945
<a href="#">gdcmBasicOffsetTable.h</a>	946
<a href="#">gdcmBitmap.h</a>	948
<a href="#">gdcmBitmapToBitmapFilter.h</a>	949
<a href="#">gdcmBoxRegion.h</a>	949
<a href="#">gdcmByteBuffer.h</a>	950
<a href="#">gdcmByteSwap.h</a>	951
<a href="#">gdcmByteSwapFilter.h</a>	952
<a href="#">gdcmByteValue.h</a>	953
<a href="#">gdcmCAPICryptoFactory.h</a>	954
<a href="#">gdcmCAPICryptographicMessageSyntax.h</a>	954
<a href="#">gdcmCEchoMessages.h</a>	955
<a href="#">gdcmCFindMessages.h</a>	956
<a href="#">gdcmCMoveMessages.h</a>	957
<a href="#">gdcmCodec.h</a>	957

<a href="#">gdcmCoder.h</a>	958
<a href="#">gdcmCodeString.h</a>	959
<a href="#">gdcmCommand.h</a>	960
<a href="#">gdcmCommandDataSet.h</a>	962
<a href="#">gdcmCompositeMessageFactory.h</a>	962
<a href="#">gdcmCompositeNetworkFunctions.h</a>	963
<a href="#">gdcmConstCharWrapper.h</a>	964
<a href="#">gdcmCP246ExplicitDataElement.h</a>	965
<a href="#">gdcmCryptoFactory.h</a>	965
<a href="#">gdcmCryptographicMessageSyntax.h</a>	966
<a href="#">gdcmCSAElement.h</a>	967
<a href="#">gdcmCSAHeader.h</a>	968
<a href="#">gdcmCSAHeaderDict.h</a>	969
<a href="#">gdcmCSAHeaderDictEntry.h</a>	971
<a href="#">gdcmCStoreMessages.h</a>	972
<a href="#">gdcmCurve.h</a>	973
<a href="#">gdcmDataElement.h</a>	974
<a href="#">gdcmDataEvent.h</a>	975
<a href="#">gdcmDataSet.h</a>	976
<a href="#">gdcmDataSetEvent.h</a>	977
<a href="#">gdcmDataSetHelper.h</a>	978
<a href="#">gdcmDecoder.h</a>	979
<a href="#">gdcmDefinedTerms.h</a>	980
<a href="#">gdcmDeflateStream.h</a>	981
<a href="#">gdcmDefs.h</a>	981
<a href="#">gdcmDeltaEncodingCodec.h</a>	983
<a href="#">gdcmDICOMDIR.h</a>	983
<a href="#">gdcmDICOMDIRGenerator.h</a>	984
<a href="#">gdcmDict.h</a>	985
<a href="#">gdcmDictConverter.h</a>	987
<a href="#">gdcmDictEntry.h</a>	987
<a href="#">gdcmDictPrinter.h</a>	989
<a href="#">gdcmDicts.h</a>	989
<a href="#">gdcmDIMSE.h</a>	991
<a href="#">gdcmDirectionCosines.h</a>	991
<a href="#">gdcmDirectory.h</a>	992
<a href="#">gdcmDirectoryHelper.h</a>	993
<a href="#">gdcmDummyValueGenerator.h</a>	994
<a href="#">gdcmDumper.h</a>	994
<a href="#">gdcmElement.h</a>	995
<a href="#">gdcmEncapsulatedDocument.h</a>	997
<a href="#">gdcmEnumeratedValues.h</a>	997
<a href="#">gdcmEvent.h</a>	998
<a href="#">gdcmException.h</a>	1000
<a href="#">gdcmExplicitDataElement.h</a>	1000
<a href="#">gdcmExplicitImplicitDataElement.h</a>	1001
<a href="#">gdcmFiducials.h</a>	1002
<a href="#">gdcmFile.h</a>	1003
<a href="#">gdcmFileAnonymizer.h</a>	1004
<a href="#">gdcmFileChangeTransferSyntax.h</a>	1004
<a href="#">gdcmFileDerivation.h</a>	1005
<a href="#">gdcmFileExplicitFilter.h</a>	1006
<a href="#">gdcmFileMetaInformation.h</a>	1007
<a href="#">gdcmFilename.h</a>	1008

gdcmFileNameEvent.h	1009
gdcmFilenameGenerator.h	1010
gdcmFileSet.h	1010
gdcmFileStreamer.h	1012
gdcmFindPatientRootQuery.h	1012
gdcmFindStudyRootQuery.h	1013
gdcmFragment.h	1014
gdcmGlobal.h	1016
gdcmGroupDict.h	1017
gdcmIconImage.h	1017
gdcmIconImageFilter.h	1018
gdcmIconImageGenerator.h	1019
gdcmImage.h	1020
gdcmImageApplyLookupTable.h	1021
gdcmImageChangePhotometricInterpretation.h	1021
gdcmImageChangePlanarConfiguration.h	1022
gdcmImageChangeTransferSyntax.h	1023
gdcmImageCodec.h	1024
gdcmImageConverter.h	1025
gdcmImageFragmentSplitter.h	1026
gdcmImageHelper.h	1027
gdcmImageReader.h	1028
gdcmImageRegionReader.h	1028
gdcmImageToImageFilter.h	1029
gdcmImageWriter.h	1030
gdcmImplementationClassUIDSub.h	1031
gdcmImplementationUIDSub.h	1032
gdcmImplementationVersionNameSub.h	1033
gdcmImplicitDataElement.h	1034
gdcmIOD.h	1035
gdcmIODEntry.h	1036
gdcmIODs.h	1038
gdcmIPPSorter.h	1039
gdcmItem.h	1040
gdcmJPEG12Codec.h	1042
gdcmJPEG16Codec.h	1042
gdcmJPEG2000Codec.h	1043
gdcmJPEG8Codec.h	1044
gdcmJPEGCodec.h	1045
gdcmJPEGLSCCodec.h	1046
gdcmJSON.h	1047
gdcmKAKADUCCodec.h	1048
gdcmLegacyMacro.h	1048
gdcmLO.h	1049
gdcmLookupTable.h	1050
gdcmMacro.h	1051
gdcmMacroEntry.h	1053
gdcmMacros.h	1054
gdcmMaximumLengthSub.h	1056
gdcmMD5.h	1057
gdcmMediaStorage.h	1058
gdcmMeshPrimitive.h	1059
gdcmModule.h	1060
gdcmModuleEntry.h	1062

<a href="#">gdcmModules.h</a>	1063
<a href="#">gdcmMovePatientRootQuery.h</a>	1065
<a href="#">gdcmMoveStudyRootQuery.h</a>	1066
<a href="#">gdcmNestedModuleEntries.h</a>	1066
<a href="#">gdcmNetworkEvents.h</a>	1068
<a href="#">gdcmNetworkStateID.h</a>	1069
<a href="#">gdcmObject.h</a>	1070
<a href="#">gdcmOpenSSLCryptoFactory.h</a>	1070
<a href="#">gdcmOpenSSLCryptographicMessageSyntax.h</a>	1071
<a href="#">gdcmOpenSSLP7CryptoFactory.h</a>	1072
<a href="#">gdcmOpenSSLP7CryptographicMessageSyntax.h</a>	1073
<a href="#">gdcmOrientation.h</a>	1075
<a href="#">gdcmOverlay.h</a>	1075
<a href="#">gdcmParseException.h</a>	1076
<a href="#">gdcmParser.h</a>	1078
<a href="#">gdcmPatient.h</a>	1078
<a href="#">gdcmPDataTFPDU.h</a>	1079
<a href="#">gdcmPDBelement.h</a>	1080
<a href="#">gdcmPDBHeader.h</a>	1082
<a href="#">gdcmPDFCodec.h</a>	1082
<a href="#">gdcmPDUFactory.h</a>	1083
<a href="#">gdcmPersonName.h</a>	1084
<a href="#">gdcmPGXCodec.h</a>	1084
<a href="#">gdcmPhotometricInterpretation.h</a>	1085
<a href="#">gdcmPixelFormat.h</a>	1086
<a href="#">gdcmPixmap.h</a>	1087
<a href="#">gdcmPixmapReader.h</a>	1088
<a href="#">gdcmPixmapToPixmapFilter.h</a>	1089
<a href="#">gdcmPixmapWriter.h</a>	1090
<a href="#">gdcmPNMCodec.h</a>	1091
<a href="#">gdcmPreamble.h</a>	1092
<a href="#">gdcmPresentationContext.h</a>	1093
<a href="#">gdcmPresentationContextAC.h</a>	1094
<a href="#">gdcmPresentationContextGenerator.h</a>	1095
<a href="#">gdcmPresentationContextRQ.h</a>	1096
<a href="#">gdcmPresentationDataValue.h</a>	1097
<a href="#">gdcmPrinter.h</a>	1097
<a href="#">gdcmPrivateTag.h</a>	1098
<a href="#">gdcmProgressEvent.h</a>	1100
<a href="#">gdcmPVRGCodec.h</a>	1100
<a href="#">gdcmPythonFilter.h</a>	1101
<a href="#">gdcmQueryBase.h</a>	1102
<a href="#">gdcmQueryFactory.h</a>	1103
<a href="#">gdcmQueryImage.h</a>	1104
<a href="#">gdcmQueryPatient.h</a>	1105
<a href="#">gdcmQuerySeries.h</a>	1106
<a href="#">gdcmQueryStudy.h</a>	1107
<a href="#">gdcmRAWCodec.h</a>	1107
<a href="#">gdcmReader.h</a>	1108
<a href="#">gdcmRegion.h</a>	1110
<a href="#">gdcmRescaler.h</a>	1111
<a href="#">gdcmRLECodec.h</a>	1111
<a href="#">gdcmRoleSelectionSub.h</a>	1112
<a href="#">gdcmScanner.h</a>	1113

gdcmSegment.h	1114
gdcmSegmentedPaletteColorLookupTable.h	1115
gdcmSegmentHelper.h	1116
gdcmSegmentReader.h	1117
gdcmSegmentWriter.h	1118
gdcmSequenceOfFragments.h	1119
gdcmSequenceOfItems.h	1120
gdcmSerieHelper.h	1121
gdcmSeries.h	1122
gdcmServiceClassApplicationInformation.h	1124
gdcmServiceClassUser.h	1125
gdcmSHA1.h	1125
gdcmSimpleSubjectWatcher.h	1126
gdcmSmartPointer.h	1127
gdcmSOPClassExtendedNegociationSub.h	1128
gdcmSOPClassUIDToIOD.h	1129
gdcmSorter.h	1130
gdcmSpacing.h	1132
gdcmSpectroscopy.h	1132
gdcmSplitMosaicFilter.h	1133
gdcmStaticAssert.h	1134
gdcmStreamImageReader.h	1135
gdcmStreamImageWriter.h	1135
gdcmString.h	1136
gdcmStringFilter.h	1137
gdcmStudy.h	1138
gdcmSubject.h	1139
gdcmSurface.h	1139
gdcmSurfaceHelper.h	1141
gdcmSurfaceReader.h	1141
gdcmSurfaceWriter.h	1142
gdcmSwapCode.h	1143
gdcmSwapper.h	1144
gdcmSystem.h	1145
gdcmTable.h	1146
gdcmTableEntry.h	1146
gdcmTableReader.h	1147
gdcmTag.h	1149
gdcmTagPath.h	1149
gdcmTagToVR.h	1150
gdcmTerminal.h	1151
gdcmTestDriver.h	1153
gdcmTesting.h	1153
gdcmTrace.h	1154
gdcmTransferSyntax.h	1157
gdcmTransferSyntaxSub.h	1158
gdcmType.h	1159
gdcmTypes.h	1160
gdcmUIDGenerator.h	1161
gdcmUIDs.h	1161
gdcmULAction.h	1162
gdcmULActionAA.h	1163
gdcmULActionAE.h	1164
gdcmULActionAR.h	1165

gdcmULActionDT.h	1166
gdcmULBasicCallback.h	1166
gdcmULConnection.h	1167
gdcmULConnectionCallback.h	1168
gdcmULConnectionInfo.h	1169
gdcmULConnectionManager.h	1171
gdcmULEvent.h	1171
gdcmULTransitionTable.h	1172
gdcmULWritingCallback.h	1174
gdcmUNExplicitDataElement.h	1174
gdcmUNExplicitImplicitDataElement.h	1175
gdcmUnpacker12Bits.h	1175
gdcmUsage.h	1176
gdcmUserInformation.h	1178
gdcmUUIDGenerator.h	1179
gdcmValidate.h	1179
gdcmValue.h	1180
gdcmValueIO.h	1181
gdcmVersion.h	1182
gdcmVL.h	1183
gdcmVM.h	1184
gdcmVR.h	1185
gdcmVR16ExplicitDataElement.h	1187
gdcmWaveform.h	1188
gdcmWin32.h	1188
gdcmWriter.h	1189
gdcmXMLDictReader.h	1190
gdcmXMLPrinter.h	1190
gdcmXMLPrivateDictReader.h	1191
vtkGDCMImageReader.h	1192
vtkGDCMImageReader2.h	1193
vtkGDCMImageWriter.h	1194
vtkGDCMMedicalImageProperties.h	1194
vtkGDCMPolyDataReader.h	1195
vtkGDCMPolyDataWriter.h	1196
vtkGDCMTesting.h	1196
vtkGDCMThreadedImageReader.h	1197
vtkGDCMThreadedImageReader2.h	1198
vtkImageColorViewer.h	1198
vtkImageMapToColors16.h	1199
vtkImageMapToWindowLevelColors2.h	1199
vtkImagePlanarComponentsToComponents.h	1200
vtkImageRGBToYBR.h	1200
vtkImageYBRToRGB.h	1201
vtkLookupTable16.h	1201
vtkRTStructSetProperties.h	1202



## Chapter 26

# Namespace Documentation

### 26.1 gdcmm Namespace Reference

#### Namespaces

- [network](#)
- [SegmentHelper](#)
- [terminal](#)

*Class for Terminal Allow one to print in color in a shell.*

#### Classes

- class [AbortEvent](#)
- class [AnonymizeEvent](#)  
*[AnonymizeEvent](#) Special type of event triggered during the Anonymization process.*
- class [Anonymizer](#)  
*[Anonymizer](#) This class is a multi purpose anonymizer. It can work in 2 mode:*
- class [AnyEvent](#)
- class [ApplicationEntity](#)  
*[ApplicationEntity](#).*
- class [ASN1](#)  
*Class for [ASN1](#).*
- class [Attribute](#)  
*[Attribute](#) class This class use template metaprograming tricks to let the user know when the template instanciation does not match the public dictionary.*
- class [Attribute< Group, Element, TVR, VM::VM1 >](#)
- class [Attribute< Group, Element, TVR, VM::VM1\\_3 >](#)
- class [Attribute< Group, Element, TVR, VM::VM1\\_8 >](#)
- class [Attribute< Group, Element, TVR, VM::VM1\\_n >](#)
- class [Attribute< Group, Element, TVR, VM::VM2\\_2n >](#)
- class [Attribute< Group, Element, TVR, VM::VM2\\_n >](#)
- class [Attribute< Group, Element, TVR, VM::VM3\\_3n >](#)
- class [Attribute< Group, Element, TVR, VM::VM3\\_n >](#)
- class [AudioCodec](#)

- [AudioCodec](#).
- class [Base64](#)
  - Class for [Base64](#).*
- class [BaseRootQuery](#)
  - [BaseRootQuery](#) contains: a baseclass which will produce a dataset for c-find and c-move with patient/study root.*
- class [BasicOffsetTable](#)
  - Class to represent a [BasicOffsetTable](#).*
- class [Bitmap](#)
  - [Bitmap](#) class A bitmap based image. Used as parent for both [IconImage](#) and the main Pixel Data [Image](#) It does not contains any World Space information (IPP, IOP)*
- class [BitmapToBitmapFilter](#)
  - [BitmapToBitmapFilter](#) class Super class for all filter taking an image and producing an output image.*
- class [BoxRegion](#)
  - Class for manipulation box region This is a very simple implementation of the [Region](#) class. It only support 3D box type region. It assumes the 3D Box does not have a tilt Origin is as (0,0,0)*
- class [ByteBuffer](#)
  - [ByteBuffer](#).*
- class [ByteSwap](#)
  - [ByteSwap](#).*
- class [ByteSwapFilter](#)
  - [ByteSwapFilter](#) In place byte-swapping of a dataset FIXME: FL status ??*
- class [ByteValue](#)
  - Class to represent binary value (array of bytes)*
- class [CAPICryptoFactory](#)
- class [CAPICryptographicMessageSyntax](#)
- class [Codec](#)
  - [Codec](#) class.*
- class [Coder](#)
  - [Coder](#).*
- class [CodeString](#)
  - [CodeString](#) This is an implementation of DICOM [VR](#): CS The ctor will properly Trim so that operator== is correct.*
- class [Command](#)
  - [Command](#) superclass for callback/observer methods.*
- class [CommandDataSet](#)
  - Class to represent a [Command DataSet](#).*
- class [CompositeNetworkFunctions](#)
  - Composite Network Functions These functions provide a generic API to the DICOM functions implemented in GDCM. Advanced users can use this code as a template for building their own versions of these functions (for instance, to provide progress bars or some other way of handling returned query information), but for most users, these functions should be sufficient to interface with a PACS to a local machine. Note that these functions are not contained within a static class or some other class-style interface, because multiple connections can be instantiated in the same program. The DICOM standard is much more function oriented rather than class oriented in this instance, so the design of this API reflects that functional approach. These functions implements the following SCU operations:*
- class [ConstCharWrapper](#)
  - Do not use me.*
- class [CP246ExplicitDataElement](#)
  - Class to read/write a [DataElement](#) as CP246Explicit Data [Element](#).*
- class [CryptoFactory](#)
  - Class to do handle the crypto factory.*

- class [CryptographicMessageSyntax](#)
- class [CSAElement](#)
  - Class to represent a CSA [Element](#).*
- class [CSAHeader](#)
  - Class for [CSAHeader](#).*
- class [CSAHeaderDict](#)
  - Class to represent a map of [CSAHeaderDictEntry](#).*
- class [CSAHeaderDictEntry](#)
  - Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcM::Tag](#) to the needed information.*
- class [CSAHeaderDictException](#)
- class [Curve](#)
  - [Curve](#) class to handle element 50xx,3000 [Curve](#) Data WARNING: This is deprecated and lastly defined in PS 3.3 - 2004.*
- class [DataElement](#)
  - Class to represent a Data [Element](#) either Implicit or Explicit.*
- class [DataElementException](#)
- class [DataEvent](#)
  - [DataEvent](#).*
- class [DataSet](#)
  - Class to represent a Data Set (which contains Data Elements) A Data Set represents an instance of a real world Information [Object](#).*
- class [DataSetEvent](#)
  - [DataSetEvent](#) Special type of event triggered during the [DataSet](#) store/move process.*
- class [DataSetHelper](#)
  - [DataSetHelper](#) (internal class, not intended for user level)*
- class [Decoder](#)
  - [Decoder](#).*
- class [DefinedTerms](#)
  - Defined Terms are used when the specified explicit Values may be extended by implementors to include additional new Values. These new Values shall be specified in the Conformance Statement (see PS 3.2) and shall not have the same meaning as currently defined Values in this standard. A Data [Element](#) with Defined Terms that does not contain a [Value](#) equivalent to one of the Values currently specified in this standard shall not be considered to have an invalid value. Note: Interpretation [Type](#) ID (4008,0210) is an example of a Data [Element](#) having Defined Terms. It is defined to have a [Value](#) that may be one of the set of standard Values; REPORT or AMENDMENT (see PS 3.3). Because this Data [Element](#) has Defined Terms other Interpretation [Type](#) IDs may be defined by the implementor.*
- class [Defs](#)
  - FIXME I do not like the name '[Defs](#)'.*
- class [DeltaEncodingCodec](#)
  - [DeltaEncodingCodec](#) compression used by some private vendor.*
- class [DICOMDIR](#)
  - [DICOMDIR](#) class.*
- class [DICOMDIRGenerator](#)
  - [DICOMDIRGenerator](#) class This is a STD-GEN-CD [DICOMDIR](#) generator. ref: PS 3.11-2008 Annex D (Normative) - General Purpose CD-R and DVD Interchange Profiles.*
- class [Dict](#)
  - Class to represent a map of [DictEntry](#).*
- class [DictConverter](#)
  - Class to convert a .dic file into something else:*
- class [DictEntry](#)

Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcm::Tag](#) to the needed information.

- class [DictPrinter](#)

[DictPrinter](#) class.

- class [Dicts](#)

Class to manipulate the sum of knowledge (all the dict user load)

- class [DirectionCosines](#)

class to handle [DirectionCosines](#)

- class [Directory](#)

Class for manipulation directories.

- class [DirectoryHelper](#)

[DirectoryHelper](#) this class is designed to help mitigate some of the commonly performed operations on directories. namely: 1) the ability to determine the number of series in a directory by what type of series is present 2) the ability to find all ct series in a directory 3) the ability to find all mr series in a directory 4) to load a set of DataSets from a series that's already been sorted by the IPP sorter 5) For rtstruct stuff, you need to know the sopinstanceuid of each z plane, so there's a retrieval function for that 6) then a few other functions for rtstruct writeouts.

- class [DummyValueGenerator](#)

Class for generating dummy value.

- class [Dumper](#)

[Codec](#) class.

- class [Element](#)

[Element](#) class.

- class [Element< TVR, VM::VM1\\_2 >](#)
- class [Element< TVR, VM::VM1\\_n >](#)
- class [Element< TVR, VM::VM2\\_2n >](#)
- class [Element< TVR, VM::VM2\\_n >](#)
- class [Element< TVR, VM::VM3\\_3n >](#)
- class [Element< TVR, VM::VM3\\_n >](#)
- class [Element< VR::AS, VM::VM5 >](#)
- class [Element< VR::OB, VM::VM1 >](#)
- class [Element< VR::OW, VM::VM1 >](#)
- class [ElementDisableCombinations](#)

A class which is used to produce compile errors for an invalid combination of template parameters.

- class [ElementDisableCombinations< VR::OB, VM::VM1\\_n >](#)
- class [ElementDisableCombinations< VR::OW, VM::VM1\\_n >](#)
- class [EncapsulatedDocument](#)

[EncapsulatedDocument](#).

- class [EncodingImplementation](#)

[EncodingImplementation](#).

- class [EncodingImplementation< VR::VRASCII >](#)
- class [EncodingImplementation< VR::VRBINARY >](#)
- class [EndEvent](#)
- class [EnumeratedValues](#)

[Element](#). A Data [Element](#) with Enumerated Values that does not have a [Value](#) equivalent to one of the Values specified in this standard has an invalid value within the scope of a specific Information Object/SOP Class definition. Note:

- class [Event](#)

superclass for callback/observer methods

- class [Exception](#)

[Exception](#).

- class [ExitEvent](#)
- class [ExplicitDataElement](#)  
*Class to read/write a [DataElement](#) as Explicit Data [Element](#).*
- class [ExplicitImplicitDataElement](#)  
*Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#).*
- class [Fiducials](#)  
*Fiducials.*
- class [File](#)  
*a DICOM [File](#) See PS 3.10 [File](#): A [File](#) is an ordered string of zero or more bytes, where the first byte is at the beginning of the file and the last byte at the end of the [File](#). Files are identified by a unique [File](#) ID and may be written, read and/or deleted.*
- class [FileAnonymizer](#)  
*FileAnonymizer.*
- class [FileChangeTransferSyntax](#)  
*FileChangeTransferSyntax.*
- class [FileDerivation](#)  
*FileDerivation class See PS 3.16 - 2008 For the list of Code [Value](#) that can be used for in Derivation Code Sequence.*
- class [FileExplicitFilter](#)  
*FileExplicitFilter class After changing a file from Implicit to Explicit representation (see [ImageChangeTransferSyntax](#)) one operation is to make sure the [VR](#) of each DICOM attribute are accurate and do match the one from PS 3.6. Indeed when a file is written in Implicit representation, the [VR](#) is not stored directly in the file.*
- class [FileMetaInformation](#)  
*Class to represent a [File](#) Meta Information.*
- class [Filename](#)  
*Class to manipulate file name's.*
- class [FileNameEvent](#)  
*FileNameEvent Special type of event triggered during processing of [FileSet](#).*
- class [FilenameGenerator](#)  
*FilenameGenerator.*
- class [FileSet](#)  
*File-set: A File-set is a collection of DICOM Files (and possibly non-DICOM Files) that share a common naming space within which [File](#) IDs are unique.*
- class [FileStreamer](#)  
*FileStreamer This class let a user create a massive DICOM [DataSet](#) from a template DICOM file, by appending chunks of data.*
- class [FileWithName](#)  
*FileWithName.*
- class [FindPatientRootQuery](#)  
*PatientRootQuery contains: the class which will produce a dataset for c-find with patient root.*
- class [FindStudyRootQuery](#)  
*FindStudyRootQuery contains: the class which will produce a dataset for C-FIND with study root.*
- class [Fragment](#)  
*Class to represent a [Fragment](#).*
- class [Global](#)  
*Global.*
- class [GroupDict](#)  
*Class to represent the mapping from group number to its abbreviation and name.*
- class [IconImageFilter](#)

*IconImageFilter* This filter will extract icons from a *File* This filter will loop over all known sequence (public and private) that may contains an *IconImage* and retrieve them. The filter will fails with a value of false if no icon can be found Since it handle both public and private icon type, one should not assume the icon is in uncompress form, some private vendor store private icon in JPEG8/JPEG12.

- class *IconImageGenerator*

*IconImageGenerator* This filter will generate a valid Icon from the Pixel Data element (an instance of *Pixmap*). To generate a valid Icon, one is only allowed the following Photometric Interpretation:

- struct *ignore\_char*

- class *Image*

*Image* This is the container for an *Image* in the general sense. From this container you should be able to request information like:

- class *ImageApplyLookupTable*

*ImageApplyLookupTable* class It applies the LUT the PixelData (only PALETTE\_COLOR images) Output will be a *PhotometricInterpretation=RGB* image.

- class *ImageChangePhotometricInterpretation*

*ImageChangePhotometricInterpretation* class Class to change the Photometric Interpretation of an input DICOM.

- class *ImageChangePlanarConfiguration*

*ImageChangePlanarConfiguration* class Class to change the Planar configuration of an input DICOM By default it will change into the more usual representation: *PlanarConfiguration = 0*.

- class *ImageChangeTransferSyntax*

*ImageChangeTransferSyntax* class Class to change the transfer syntax of an input DICOM.

- class *ImageCodec*

*ImageCodec*.

- class *ImageConverter*

*Image* Converter.

- class *ImageFragmentSplitter*

*ImageFragmentSplitter* class For single frame image, DICOM standard allow splitting the frame into multiple fragments.

- class *ImageHelper*

*ImageHelper* (internal class, not intended for user level)

- class *ImageReader*

*ImageReader*.

- class *ImageRegionReader*

*ImageRegionReader*.

- class *ImageToImageFilter*

*ImageToImageFilter* class Super class for all filter taking an image and producing an output image.

- class *ImageWriter*

*ImageWriter*.

- class *ImplicitDataElement*

Class to represent an Implicit VR Data Element.

- class *InitializeEvent*

- class *IOD*

Class for representing a IOD.

- class *IODEntry*

Class for representing a IODEntry.

- class *IODs*

Class for representing a IODs.

- class *IPPSorter*

*IPPSorter* Implement a simple *Image* Position (*Patient*) sorter, along the *Image Orientation* (*Patient*) direction. This algorithm does NOT support duplicate and will FAIL in case of duplicate IPP.

- class [Item](#)  
*Class to represent an [Item](#) A component of the value of a Data [Element](#) that is of [Value](#) Representation Sequence of Items. An [Item](#) contains a Data Set . See PS 3.5 7.5.1 [Item](#) Encoding Rules Each [Item](#) of a Data [Element](#) of VR SQ shall be encoded as a DICOM Standart Data [Element](#) with a specific Data [Element](#) Tag of [Value](#) (FFFE,E000). The [Item](#) Tag is followed by a 4 byte [Item](#) Length field encoded in one of the following two ways Explicit/ Implicit.*
- class [IterationEvent](#)
- class [JPEG12Codec](#)  
*Class to do JPEG 12bits (lossy & lossless)*
- class [JPEG16Codec](#)  
*Class to do JPEG 16bits (lossless)*
- class [JPEG2000Codec](#)  
*Class to do JPEG 2000.*
- class [JPEG8Codec](#)  
*Class to do JPEG 8bits (lossy & lossless)*
- class [JPEGCodec](#)  
*JPEG codec Class to do JPEG (8bits, 12bits, 16bits lossy & lossless). It redispach in between the different codec implementation: [JPEG8Codec](#), [JPEG12Codec](#) & [JPEG16Codec](#) It also support inconsistency in between DICOM header and JPEG compressed stream [ImageCodec](#) implementation for the JPEG case.*
- class [JPEGLSCCodec](#)  
*JPEG-LS.*
- class [JSON](#)
- class [KAKADUCodec](#)  
*KAKADUCodec.*
- class [LO](#)  
*LO.*
- class [LookupTable](#)  
*[LookupTable](#) class.*
- class [Macro](#)  
*Class for representing a [Macro](#).*
- class [Macros](#)  
*Class for representing a [Modules](#).*
- class [MD5](#)  
*Class for [MD5](#).*
- class [MediaStorage](#)  
*[MediaStorage](#).*
- class [MemberCommand](#)  
*[Command](#) subclass that calls a pointer to a member function.*
- class [MeshPrimitive](#)  
*This class defines surface mesh primitives. It is designed from surface mesh primitives macro.*
- class [ModifiedEvent](#)
- class [Module](#)  
*Class for representing a [Module](#).*
- class [ModuleEntry](#)  
*Class for representing a [ModuleEntry](#).*
- class [Modules](#)  
*Class for representing a [Modules](#).*
- class [MovePatientRootQuery](#)  
*[MovePatientRootQuery](#) contains: the class which will produce a dataset for c-move with patient root.*

- class [MoveStudyRootQuery](#)  
*MoveStudyRootQuery* contains: the class which will produce a dataset for C-MOVE with study root.
- class [NestedModuleEntries](#)  
*Class for representing a [NestedModuleEntries](#).*
- class [NoEvent](#)
- class [Object](#)  
*Object.*
- class [OpenSSLCryptoFactory](#)
- class [OpenSSLCryptographicMessageSyntax](#)
- class [OpenSSL7CryptoFactory](#)
- class [OpenSSL7CryptographicMessageSyntax](#)  
*Class for [CryptographicMessageSyntax](#) encryption. This is just a simple wrapper around openssl PKCS7\_encrypt functionalities.*
- class [Orientation](#)  
*class to handle [Orientation](#)*
- class [Overlay](#)  
*Overlay class.*
- class [ParseException](#)  
*ParseException* Standard exception handling object.
- class [Parser](#)  
*Parser* ala XML\_Parser from expat (SAX)
- class [Patient](#)  
*See PS 3.3 - 2007 DICOM MODEL OF THE REAL-WORLD, p 54.*
- class [PDBElement](#)  
*Class to represent a PDB [Element](#).*
- class [PDBHeader](#)  
*Class for [PDBHeader](#).*
- class [PDFCodec](#)  
*PDFCodec class.*
- class [PersonName](#)  
*PersonName class.*
- class [PGXCodec](#)  
*Class to do PGX See PGX as used in JPEG 2000 implementation and reference images.*
- class [PhotometricInterpretation](#)  
*Class to represent an [PhotometricInterpretation](#).*
- class [PixelFormat](#)  
*PixelFormat.*
- class [Pixmap](#)  
*Pixmap class A bitmap based image. Used as parent for both IconImage and the main Pixel Data [Image](#) It does not contains any World Space information (IPP, IOP)*
- class [PixmapReader](#)  
*PixmapReader.*
- class [PixmapToPixmapFilter](#)  
*PixmapToPixmapFilter class Super class for all filter taking an image and producing an output image.*
- class [PixmapWriter](#)  
*PixmapWriter This class will takes two inputs:*
- class [PNMCodec](#)



Class to do PNM PNM is the Portable anymap file format. The main web page can be found at: <http://netpbm.sourceforge.net/>.

- class [Preamble](#)  
*DICOM Preamble (Part 10)*
- class [PresentationContext](#)  
*PresentationContext.*
- class [PresentationContextGenerator](#)  
*PresentationContextGenerator* This class is responsible for generating the proper [PresentationContext](#) that will be used in subsequent operation during a DICOM Query/Retrieve association. The step of the association is very sensible as special care need to be taken to explicitly define what instance are going to be send and how they are encoded.
- class [Printer](#)  
*Printer* class.
- class [PrivateDict](#)  
*Private Dict.*
- class [PrivateTag](#)  
*Class to represent a Private DICOM Data Element (Attribute) Tag (Group, Element, Owner)*
- class [ProgressEvent](#)  
*ProgressEvent* Special type of event triggered during.
- class [PVRGCodec](#)  
*PVRGCodec.*
- class [PythonFilter](#)  
*PythonFilter* *PythonFilter* is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a *DataElement* into a string, typically this is a nice feature to have for wrapped language.
- class [QueryBase](#)  
*QueryBase* contains: the base class for constructing a query dataset for a C-FIND and a C-MOVE.
- class [QueryFactory](#)  
*QueryFactory.h.*
- class [QueryImage](#)  
*QueryImage* contains: class to construct an image-based query for C-FIND and C-MOVE.
- class [QueryPatient](#)  
*QueryPatient* contains: class to construct a patient-based query for c-find and c-move.
- class [QuerySeries](#)  
*QuerySeries* contains: class to construct a series-based query for c-find and c-move.
- class [QueryStudy](#)  
*QueryStudy.h* contains: class to construct a study-based query for C-FIND and C-MOVE.
- class [RAWCodec](#)  
*RAWCodec* class.
- class [Reader](#)  
*Reader* ala DOM (Document Object Model)
- class [Region](#)  
*Class for manipulation region.*
- class [Rescaler](#)  
*Rescale class* This class is meant to apply the linear transform of Stored Pixel *Value* to Real World *Value*. This is mostly found in CT or PET dataset, where the value are stored using one type, but need to be converted to another scale using a linear transform. There are basically two cases: In CT: the linear transform is generally integer based. E.g. the Stored Pixel *Type* is unsigned short 12bits, but to get Hounsfield unit, one need to apply the linear transform:

$$RWV = 1. * SV - 1024$$

So the best scalar to store the Real World *Value* will be 16 bits signed type.

- class [RLECodec](#)  
*Class to do RLE.*
- class [Scanner](#)  
*Scanner This filter is meant for quickly browsing a [FileSet](#) (a set of files on disk). Special consideration are taken so as to read the minimum amount of information in each file in order to retrieve the user specified set of DICOM [Attribute](#).*
- class [Segment](#)  
*This class defines a segment. It mainly contains attributes of group 0x0062. In addition, it can be associated with surface.*
- class [SegmentedPaletteColorLookupTable](#)  
*SegmentedPaletteColorLookupTable class.*
- class [SegmentReader](#)  
*This class defines a segment reader. It reads attributes of group 0x0062.*
- class [SegmentWriter](#)  
*This class defines a segment writer. It writes attributes of group 0x0062.*
- class [SequenceOfFragments](#)  
*Class to represent a Sequence Of Fragments.*
- class [SequenceOfItems](#)  
*Class to represent a Sequence Of Items (value representation : SQ)*
- class [SerieHelper](#)  
*SerieHelper DO NOT USE this class, it is only a temporary solution for ITK migration from GDCM 1.x to GDCM 2.x It will disappear soon, you've been warned.*
- class [Series](#)  
*Series.*
- class [ServiceClassUser](#)  
*ServiceClassUser.*
- class [SHA1](#)  
*Class for SHA1.*
- class [SimpleMemberCommand](#)  
*Command subclass that calls a pointer to a member function.*
- class [SimpleSubjectWatcher](#)  
*SimpleSubjectWatcher This is a typical [Subject](#) Watcher class. It will observe all events.*
- class [SmartPointer](#)  
*Class for Smart Pointer.*
- class [SOPClassUIDToIOD](#)  
*Class convert a class SOP Class UID into [IOD](#).*
- class [Sorter](#)  
*Sorter General class to do sorting using a custom function You simply need to provide a function of type: [Sorter::SortFunction](#).*
- class [Spacing](#)  
*Class for Spacing.*
- class [Spectroscopy](#)  
*Spectroscopy class.*
- class [SplitMosaicFilter](#)  
*SplitMosaicFilter class Class to reshuffle bytes for a SIEMENS Mosaic image Siemens CSA [Image](#) Header CSA:= Common Siemens Architecture, sometimes also known as Common syngo Architecture.*
- class [StartEvent](#)
- struct [static\\_assert\\_test](#)
- struct [STATIC\\_ASSERTION\\_FAILURE](#)
- struct [STATIC\\_ASSERTION\\_FAILURE< true >](#)

- class [StreamImageReader](#)  
*StreamImageReader.*
- class [StreamImageWriter](#)  
*StreamImageReader.*
- class [String](#)  
*String.*
- class [StringFilter](#)  
*StringFilter* *StringFilter* is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a *DataElement* into a string, typically this is a nice feature to have for wrapped language.
- class [Study](#)  
*Study.*
- class [Subject](#)  
*Subject.*
- class [Surface](#)  
*This class defines a SURFACE IE. This members are taken from required surface mesh module attributes.*
- class [SurfaceHelper](#)  
*SurfaceHelper* Helper class for *Surface* object.
- class [SurfaceReader](#)  
*This class defines a SURFACE IE reader. It reads surface mesh module attributes.*
- class [SurfaceWriter](#)  
*This class defines a SURFACE IE writer. It writes surface mesh module attributes.*
- class [SwapCode](#)  
*SwapCode* representation.
- class [SwapperDoOp](#)
- class [SwapperNoOp](#)
- class [System](#)  
*Class to do system operation.*
- class [Table](#)  
*Table.*
- class [TableEntry](#)  
*TableEntry.*
- class [TableReader](#)  
*Class for representing a [TableReader](#).*
- class [Tag](#)  
*Class to represent a DICOM Data [Element](#) ([Attribute](#)) [Tag](#) (Group, [Element](#)). Basically an uint32\_t which can also be expressed as two uint16\_t (group and element)*
- class [TagPath](#)  
*class to handle a path of tag.*
- class [Testing](#)  
*class for testing*
- class [Trace](#)  
*Trace.*
- class [TransferSyntax](#)  
*Class to manipulate Transfer Syntax.*
- class [Type](#)  
*Type.*
- struct [UI](#)
- class [UIDGenerator](#)

- Class for generating unique UID.*
- class [UIDs](#)
  - all known uids*
- class [UNExplicitDataElement](#)
  - Class to read/write a [DataElement](#) as UNExplicit Data [Element](#).*
- class [UNExplicitImplicitDataElement](#)
  - Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#) This class gather two known bugs:*
- class [Unpacker12Bits](#)
  - Pack/Unpack 12 bits pixel into 16bits.*
- class [Usage](#)
  - Usage.*
- class [UserEvent](#)
- class [UUIDGenerator](#)
  - Class for generating unique UUID generate DCE 1.1 uid.*
- class [Validate](#)
  - [Validate](#) class.*
- class [Value](#)
  - Class to represent the value of a Data [Element](#).*
- class [ValueIO](#)
  - Class to dispatch template calls.*
- class [Version](#)
  - major/minor and build version*
- class [VL](#)
  - [Value](#) Length.*
- class [VM](#)
  - [Value](#) Multiplicity Looking at the DICOMV3 dict only there is very few cases: 1 2 3 4 5 6 8 16 24 1-2 1-3 1-8 1-32 1-99 1-n 2-2n 2-n 3-3n 3-n.*
- struct [VMToLength](#)
- class [VR](#)
  - [VR](#) class This is adapted from DICOM standard The biggest difference is the INVALID [VR](#) and the composite one that differ from standard (more like an addition) This allow us to represent all the possible case express in the DICOMV3 dict.*
- class [VR16ExplicitDataElement](#)
  - Class to read/write a [DataElement](#) as Explicit Data [Element](#).*
- struct [VRToEncoding](#)
- struct [VRToType](#)
- class [VRVLSize](#)
- class [VRVLSize< 0 >](#)
- class [VRVLSize< 1 >](#)
- class [Waveform](#)
  - [Waveform](#) class.*
- class [Writer](#)
  - [Writer](#) ala DOM (Document [Object](#) Model) This class is a non-validating writer, it will only performs well- formedness check only.*
- class [XMLDictReader](#)
  - Class for representing a [XMLDictReader](#).*
- class [XMLPrinter](#)
- class [XMLPrivateDictReader](#)
  - Class for representing a [XMLPrivateDictReader](#).*

## Typedefs

- typedef [String](#)<"\", 16 > [AECComp](#)
- typedef [String](#)<"\", 64 > [ASComp](#)
- typedef bool(\* [BOOL\\_FUNCTION\\_PFILE\\_PFILE\\_POINTER](#)) ([File](#) \*, [File](#) \*)
- typedef [String](#)<"\", 16 > [CSCComp](#)
- typedef [String](#)<"\", 64 > [DACComp](#)
- typedef [String](#)<"\", 64 > [DTComp](#)
- typedef std::vector< [SmartPointer](#)< [FileWithName](#) > > [FileList](#)
- typedef [Bitmap](#) [IconImage](#)
- typedef [String](#)<"\", 64 > [LOComp](#)
- typedef [String](#)<"\", 64 > [LTComp](#)
- typedef [ModuleEntry](#) [MacroEntry](#)
- typedef [NestedModuleEntries](#) [NestedMacroEntries](#)
- typedef [String](#)<"\", 64 > [PNComp](#)
- typedef [String](#)<"\", 64 > [SHComp](#)
- typedef [String](#)<"\", 64 > [STComp](#)
- typedef [String](#)<"\", 16 > [TMComp](#)
- typedef [String](#)<"\", 64, 0 > [UIComp](#)
- typedef [String](#)<"\", 64 > [UTComp](#)

## Enumerations

- enum [CompOperators](#) {  
[GDCM\\_EQUAL](#) = 0,  
[GDCM\\_DIFFERENT](#),  
[GDCM\\_GREATER](#),  
[GDCM\\_GREATEROREQUAL](#),  
[GDCM\\_LESS](#),  
[GDCM\\_LESOREQUAL](#) }
- enum [ECharSet](#) {  
[eLatin1](#) = 0,  
[eLatin2](#),  
[eLatin3](#),  
[eLatin4](#),  
[eCyrillic](#),  
[eArabic](#),  
[eGreek](#),  
[eHebrew](#),  
[eLatin5](#),  
[eJapanese](#),  
[eThai](#),  
[eJapaneseKanjiMultibyte](#),  
[eJapaneseSupplementaryKanjiMultibyte](#),  
[eKoreanHangulHanjaMultibyte](#),  
[eUTF8](#),  
[eGB18030](#) }
- enum [EQueryLevel](#) {  
[ePatient](#) = 0,  
[eStudy](#) = 1,  
[eSeries](#) = 2,  
[eImage](#) = 3 }

- enum [EQueryType](#) {  
  [eFind](#) = 0,  
  [eMove](#) }
- enum [ERootType](#) {  
  [ePatientRootType](#),  
  [eStudyRootType](#) }
- enum [LodModeType](#) {  
  [LD\\_ALL](#) = 0x00000000,  
  [LD\\_NOSEQ](#) = 0x00000001,  
  [LD\\_NOSHADOW](#) = 0x00000002,  
  [LD\\_NOSHADOWSEQ](#) = 0x00000004 }

## Functions

- [ignore\\_char](#) const [backslash](#) ('\\')
  - [VR::VRType](#) [GetVRFromTag](#) ([Tag](#) const &tag)
  - bool [operator!=](#) (const [CodeString](#) &ref, const [CodeString](#) &cs)
  - bool [operator!=](#) (const [DataElement](#) &lhs, const [DataElement](#) &rhs)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Version](#) &v)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [NestedModuleEntries](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [SwapCode](#) &sc)
  - std::ostream & [operator<<](#) (std::ostream &os, const [FileSet](#) &f)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Region](#) &r)
  - std::ostream & [operator<<](#) (std::ostream &os, [Event](#) &e)
- Generic inserter operator for [Event](#) and its subclasses.*
- std::ostream & [operator<<](#) (std::ostream &os, const [PDBelement](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [CommandDataSet](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Orientation](#) &o)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [IODs](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Macros](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Modules](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Type](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [ModuleEntry](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [GroupDict](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [PrivateTag](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [IOD](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [File](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Usage](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Sorter](#) &s)
  - std::ostream & [operator<<](#) (std::ostream &os, const [CSAHeaderDictEntry](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Preamble](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [IODEntry](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Macro](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [CSAHeaderDict](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Dicts](#) &d)
  - std::ostream & [operator<<](#) (std::ostream &os, const [PDBHeader](#) &d)
  - std::ostream & [operator<<](#) (std::ostream &os, const [CodeString](#) &str)
  - std::ostream & [operator<<](#) (std::ostream &\_os, const [Module](#) &\_val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [PhotometricInterpretation](#) &val)
  - std::ostream & [operator<<](#) (std::ostream &os, const [Directory](#) &d)

- `std::ostream & operator<< (std::ostream &os, const Global &g)`
- `std::ostream & operator<< (std::ostream &os, const Object &obj)`
- `std::ostream & operator<< (std::ostream &os, const BasicOffsetTable &val)`
- `std::ostream & operator<< (std::ostream &os, const DictEntry &val)`
- `std::ostream & operator<< (std::ostream &os, const VL &val)`
- `std::ostream & operator<< (std::ostream &os, const CSAElement &val)`
- `std::ostream & operator<< (std::ostream &os, const CSAHeader &d)`
- `std::ostream & operator<< (std::ostream &_os, const TransferSyntax &ts)`
- `std::ostream & operator<< (std::ostream &os, const FileMetaInformation &val)`
- `std::ostream & operator<< (std::ostream &_os, const VM &_val)`
- `std::ostream & operator<< (std::ostream &os, const Scanner &s)`
- `std::ostream & operator<< (std::ostream &os, const Dict &val)`
- `std::ostream & operator<< (std::ostream &_os, const MediaStorage &ms)`
- `std::ostream & operator<< (std::ostream &_os, const VR &val)`
- `std::ostream & operator<< (std::ostream &os, const Fragment &val)`
- `std::ostream & operator<< (std::ostream &os, const PixelFormat &pf)`
- `std::ostream & operator<< (std::ostream &_os, const UI &_val)`
- `std::ostream & operator<< (std::ostream &os, const DataElement &val)`
- `std::ostream & operator<< (std::ostream &_os, const Tag &_val)`
- `std::ostream & operator<< (std::ostream &os, const DataSet &val)`
- `std::ostream & operator<< (std::ostream &os, const Item &val)`
- `std::ostream & operator<< (std::ostream &os, const PrivateDict &val)`
- `std::ostream & operator<< (std::ostream &_os, const UIDs &uid)`
- `bool operator== (const CodeString &ref, const CodeString &cs)`
- `template<char TDelimiter, unsigned int TMaxLength, char TPadChar>  
std::istream & operator>> (std::istream &is, String< TDelimiter, TMaxLength, TPadChar > &ms)`
- `std::istream & operator>> (std::istream &in, ignore\_char const &ic)`
- `std::istream & operator>> (std::istream &_is, Tag &_val)`
- `template<typename Float >  
std::string to\_string (Float data)`
- `TYPETOENCODING (SQ, VRBINARY, unsigned char) TYPETOENCODING(UN`

## Variables

- static [Global](#) [GlobalInstance](#)
- [VRBINARY](#)

### 26.1.1 Detailed Description

This header defines the classes for the AA Actions, Association Abort Related Actions ([Table 9-9](#) of ps 3.8-2009).

Since each class is essentially a placeholder for a function pointer, I'm breaking with having each class have its own file for the sake of brevity of the number of files.

This header defines the classes for the AE Actions, Association Establishment Related Actions ([Table 9-6](#) of ps 3.8-2009).

Since each class is essentially a placeholder for a function pointer, I'm breaking with having each class have its own file for the sake of brevity of the number of files.

This header defines the classes for the AR Actions, Association Release Related Actions ([Table 9-8](#) of ps 3.8-2009).

Since each class is essentially a placeholder for a function pointer, I'm breaking with having each class have its own file for the sake of brevity of the number of files.

This header defines the classes for the DT Actions, Data Transfer Related Actions ([Table 9-8](#) of ps 3.8-2009).

Since each class is essentially a placeholder for a function pointer, I'm breaking with having each class have its own file for the sake of brevity of the number of files.

## 26.1.2 Typedef Documentation

26.1.2.1 `typedef String<'\',16> gdcm::AECComp`

26.1.2.2 `typedef String<'\',64> gdcm::ASComp`

26.1.2.3 `typedef bool(* gdcm::BOOL_FUNCTION_PFILE_PFILE_POINTER) (File *, File *)`

26.1.2.4 `typedef String<'\',16> gdcm::CSCComp`

26.1.2.5 `typedef String<'\',64> gdcm::DACComp`

26.1.2.6 `typedef String<'\',64> gdcm::DTComp`

26.1.2.7 `typedef std::vector<SmartPointer<FileWithName> > gdcm::FileList`

26.1.2.8 `typedef Bitmap gdcm::IconImage`

26.1.2.9 `typedef String<'\',64> gdcm::LOComp`

26.1.2.10 `typedef String<'\',64> gdcm::LTComp`

26.1.2.11 `typedef ModuleEntry gdcm::MacroEntry`

26.1.2.12 `typedef NestedModuleEntries gdcm::NestedMacroEntries`

26.1.2.13 `typedef String<'\',64> gdcm::PNComp`

26.1.2.14 `typedef String<'\',64> gdcm::SHComp`

26.1.2.15 `typedef String<'\',64> gdcm::STComp`

26.1.2.16 `typedef String<'\',16> gdcm::TMComp`

26.1.2.17 `typedef String<'\',64,0> gdcm::UIComp`

26.1.2.18 `typedef String<'\',64> gdcm::UTComp`

## 26.1.3 Enumeration Type Documentation

26.1.3.1 `enum gdcm::CompOperators`

Enumerator

***GDCM\_EQUAL***



***GDCM\_DIFFERENT***  
***GDCM\_GREATER***  
***GDCM\_GREATEROREQUAL***  
***GDCM\_LESS***  
***GDCM\_LESSCOREQUAL***

#### 26.1.3.2 enum gdcm::ECharSet

The character sets enumerated in PS 3.3 2009 Annex C, section C.12.1.1.2 The resulting character set is stored in 0008,0005 The conversion to the data element is performed by the [QueryFactory](#) itself

Enumerator

***eLatin1***  
***eLatin2***  
***eLatin3***  
***eLatin4***  
***eCyrillic***  
***eArabic***  
***eGreek***  
***eHebrew***  
***eLatin5***  
***eJapanese***  
***eThai***  
***eJapaneseKanjiMultibyte***  
***eJapaneseSupplementaryKanjiMultibyte***  
***eKoreanHangulHanjaMultibyte***  
***eUTF8***  
***eGB18030***

#### 26.1.3.3 enum gdcm::EQueryLevel

Enumerator

***ePatient***  
***eStudy***  
***eSeries***  
***eImage***

#### 26.1.3.4 enum gdcm::EQueryType

Enumerator

***eFind***  
***eMove***

## 26.1.3.5 enum gdcm::ERootType

Enumerator

***ePatientRootType***

***eStudyRootType***

## 26.1.3.6 enum gdcm::LodModeType

Enumerator

***LD\_ALL***

***LD\_NOSEQ***

***LD\_NOSHADOW***

***LD\_NOSHADOWSEQ***

## 26.1.4 Function Documentation

## 26.1.4.1 ignore\_char const gdcm::backslash ( '\ ' )

Referenced by gdcm::EncodingImplementation< VR::VRASCII >::ReadComputeLength().

## 26.1.4.2 VR::VRType gdcm::GetVRFromTag ( Tag const &amp; tag )

## 26.1.4.3 bool gdcm::operator!= ( const CodeString &amp; ref, const CodeString &amp; cs ) [inline]

## 26.1.4.4 bool gdcm::operator!= ( const DataElement &amp; lhs, const DataElement &amp; rhs ) [inline]

## 26.1.4.5 std::ostream&amp; gdcm::operator&lt;&lt; ( std::ostream &amp; os, const Version &amp; v ) [inline]

References gdcm::Version::Print().

## 26.1.4.6 std::ostream&amp; gdcm::operator&lt;&lt; ( std::ostream &amp; \_os, const NestedModuleEntries &amp; \_val ) [inline]

References gdcm::ModuleEntry::DataElementType, gdcm::ModuleEntry::DescriptionField, and gdcm::ModuleEntry::Name.

## 26.1.4.7 std::ostream&amp; gdcm::operator&lt;&lt; ( std::ostream &amp; os, const SwapCode &amp; sc ) [inline]

References gdcm::SwapCode::GetSwapCodeString().

## 26.1.4.8 std::ostream&amp; gdcm::operator&lt;&lt; ( std::ostream &amp; os, const FileSet &amp; f ) [inline]

## 26.1.4.9 std::ostream&amp; gdcm::operator&lt;&lt; ( std::ostream &amp; os, const Region &amp; r ) [inline]

References gdcm::Region::Print().

26.1.4.10 `std::ostream& gdcm::operator<< ( std::ostream & os, Event & e ) [inline]`

Generic inserter operator for [Event](#) and its subclasses.

References `gdcm::Event::Print()`.

26.1.4.11 `std::ostream& gdcm::operator<< ( std::ostream & os, const PDBElement & val ) [inline]`

References `gdcm::PDBElement::NameField`, and `gdcm::PDBElement::ValueField`.

26.1.4.12 `std::ostream& gdcm::operator<< ( std::ostream & os, const CommandDataSet & val ) [inline]`

References `gdcm::DataSet::Print()`.

26.1.4.13 `std::ostream& gdcm::operator<< ( std::ostream & os, const Orientation & o ) [inline]`

References `gdcm::Orientation::Print()`.

26.1.4.14 `std::ostream& gdcm::operator<< ( std::ostream & _os, const IODs & _val ) [inline]`

26.1.4.15 `std::ostream& gdcm::operator<< ( std::ostream & _os, const Macros & _val ) [inline]`

26.1.4.16 `std::ostream& gdcm::operator<< ( std::ostream & _os, const Modules & _val ) [inline]`

26.1.4.17 `std::ostream& gdcm::operator<< ( std::ostream & _os, const Type & val ) [inline]`

References `gdcm::Type::GetTypeString()`.

26.1.4.18 `std::ostream& gdcm::operator<< ( std::ostream & _os, const ModuleEntry & _val ) [inline]`

References `gdcm::ModuleEntry::DataElementType`, `gdcm::ModuleEntry::DescriptionField`, and `gdcm::ModuleEntry::Name`.

26.1.4.19 `std::ostream& gdcm::operator<< ( std::ostream & _os, const GroupDict & _val ) [inline]`

References `gdcm::GroupDict::GetAbbreviation()`, `gdcm::GroupDict::GetName()`, and `gdcm::GroupDict::Size()`.

26.1.4.20 `std::ostream& gdcm::operator<< ( std::ostream & os, const PrivateTag & val ) [inline]`

26.1.4.21 `std::ostream& gdcm::operator<< ( std::ostream & _os, const IOD & _val ) [inline]`

26.1.4.22 `std::ostream& gdcm::operator<< ( std::ostream & os, const File & val ) [inline]`

References `gdcm::File::GetHeader()`.

26.1.4.23 `std::ostream& gdcm::operator<< ( std::ostream & _os, const Usage & val ) [inline]`

References `gdcm::Usage::GetUsageString()`.

26.1.4.24 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Sorter & s )` `[inline]`

References `gdcmm::Sorter::Print()`.

26.1.4.25 `std::ostream& gdcmm::operator<< ( std::ostream & os, const CSAHeaderDictEntry & val )` `[inline]`

26.1.4.26 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Preamble & val )` `[inline]`

26.1.4.27 `std::ostream& gdcmm::operator<< ( std::ostream & _os, const IODEntry & _val )` `[inline]`

26.1.4.28 `std::ostream& gdcmm::operator<< ( std::ostream & _os, const Macro & _val )` `[inline]`

26.1.4.29 `std::ostream& gdcmm::operator<< ( std::ostream & os, const CSAHeaderDict & val )` `[inline]`

26.1.4.30 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Dicts & d )` `[inline]`

26.1.4.31 `std::ostream& gdcmm::operator<< ( std::ostream & os, const PDBHeader & d )` `[inline]`

References `gdcmm::PDBHeader::Print()`.

26.1.4.32 `std::ostream& gdcmm::operator<< ( std::ostream & os, const CodeString & str )` `[inline]`

26.1.4.33 `std::ostream& gdcmm::operator<< ( std::ostream & _os, const Module & _val )` `[inline]`

26.1.4.34 `std::ostream& gdcmm::operator<< ( std::ostream & os, const PhotometricInterpretation & val )` `[inline]`

References `gdcmm::PhotometricInterpretation::GetPIString()`.

26.1.4.35 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Directory & d )` `[inline]`

References `gdcmm::Directory::Print()`.

26.1.4.36 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Global & g )` `[inline]`

26.1.4.37 `std::ostream& gdcmm::operator<< ( std::ostream & os, const Object & obj )` `[inline]`

References `gdcmm::Object::Print()`.

26.1.4.38 `std::ostream& gdcmm::operator<< ( std::ostream & os, const BasicOffsetTable & val )` `[inline]`

References `gdcmm::DataElement::GetByteValue()`, `gdcmm::DataElement::ValueField`, and `gdcmm::DataElement::Value↔LengthField`.

26.1.4.39 `std::ostream& gdcmm::operator<< ( std::ostream & os, const DictEntry & val )` `[inline]`

26.1.4.40 `std::ostream& gdcmm::operator<< ( std::ostream & os, const VL & val )` `[inline]`

26.1.4.41 `std::ostream& gdcm::operator<< ( std::ostream & os, const CSAElement & val ) [inline]`

References `gdcm::CSAElement::DataField`, `gdcm::ByteValue::GetLength()`, `gdcm::ByteValue::GetPointer()`, `gdcm::CSAElement::KeyField`, `gdcm::CSAElement::NameField`, `gdcm::CSAElement::NoOfItemsField`, `gdcm::CSAElement::SyngoDTField`, `gdcm::CSAElement::ValueMultiplicityField`, `gdcm::VM::VM1`, and `gdcm::CSAElement::VRField`.

26.1.4.42 `std::ostream& gdcm::operator<< ( std::ostream & os, const CSAHeader & d ) [inline]`

References `gdcm::CSAHeader::Print()`.

26.1.4.43 `std::ostream& gdcm::operator<< ( std::ostream & _os, const TransferSyntax & ts ) [inline]`

References `gdcm::TransferSyntax::GetTSSString()`.

26.1.4.44 `std::ostream& gdcm::operator<< ( std::ostream & os, const FileMetaInformation & val ) [inline]`

References `gdcm::FileMetaInformation::GetPreamble()`, and `gdcm::DataSet::Print()`.

26.1.4.45 `std::ostream& gdcm::operator<< ( std::ostream & _os, const VM & _val ) [inline]`

References `gdcm::VM::GetVMString()`.

26.1.4.46 `std::ostream& gdcm::operator<< ( std::ostream & os, const Scanner & s ) [inline]`

References `gdcm::Scanner::Print()`.

26.1.4.47 `std::ostream& gdcm::operator<< ( std::ostream & os, const Dict & val ) [inline]`

26.1.4.48 `std::ostream& gdcm::operator<< ( std::ostream & _os, const MediaStorage & ms ) [inline]`

References `gdcm::MediaStorage::GetMSString()`.

26.1.4.49 `std::ostream& gdcm::operator<< ( std::ostream & _os, const VR & val ) [inline]`

References `gdcm::VR::GetVRString()`.

26.1.4.50 `std::ostream& gdcm::operator<< ( std::ostream & os, const Fragment & val ) [inline]`

References `gdcm::DataElement::TagField`, `gdcm::DataElement::ValueField`, and `gdcm::DataElement::ValueLengthField`.

26.1.4.51 `std::ostream& gdcm::operator<< ( std::ostream & os, const PixelFormat & pf ) [inline]`

References `gdcm::PixelFormat::Print()`.

26.1.4.52 `std::ostream& gdcm::operator<< ( std::ostream & _os, const UI & _val )` [inline]

References `gdcm::UI::Internal`.

26.1.4.53 `std::ostream& gdcm::operator<< ( std::ostream & os, const DataElement & val )` [inline]

References `gdcm::Object::Print()`, `gdcm::DataElement::TagField`, `gdcm::DataElement::ValueField`, `gdcm::DataElement::ValueLengthField`, and `gdcm::DataElement::VRField`.

26.1.4.54 `std::ostream& gdcm::operator<< ( std::ostream & _os, const Tag & _val )` [inline]

26.1.4.55 `std::ostream& gdcm::operator<< ( std::ostream & os, const DataSet & val )` [inline]

References `gdcm::DataSet::Print()`.

26.1.4.56 `std::ostream& gdcm::operator<< ( std::ostream & os, const Item & val )` [inline]

References `gdcm::DataSet::Print()`, `gdcm::DataElement::TagField`, and `gdcm::DataElement::ValueLengthField`.

26.1.4.57 `std::ostream& gdcm::operator<< ( std::ostream & os, const PrivateDict & val )` [inline]

26.1.4.58 `std::ostream& gdcm::operator<< ( std::ostream & _os, const UIDs & uid )` [inline]

References `gdcm::UIDs::GetName()`, and `gdcm::UIDs::GetString()`.

26.1.4.59 `bool gdcm::operator== ( const CodeString & ref, const CodeString & cs )` [inline]

Examples:

[DumpPhilipsECHO.cxx](#).

26.1.4.60 `template<char TDelimiter, unsigned int TMaxLength, char TPadChar> std::istream& gdcm::operator>> ( std::istream & is, String< TDelimiter, TMaxLength, TPadChar > & ms )` [inline]

26.1.4.61 `std::istream& gdcm::operator>> ( std::istream & in, ignore_char const & ic )` [inline]

References `gdcm::ignore_char::m_char`.

26.1.4.62 `std::istream& gdcm::operator>> ( std::istream & _is, Tag & _val )` [inline]

References `gdcm::Tag::SetElement()`, and `gdcm::Tag::SetGroup()`.

26.1.4.63 `template<typename Float > std::string gdcm::to_string ( Float data )`

Referenced by `gdcm::EncodingImplementation< VR::VRASCII >::Write()`.

26.1.4.64 gdcmm::TYPETOENCODING ( SQ , VRBINARY , unsigned char )

## 26.1.5 Variable Documentation

26.1.5.1 Global gdcmm::GlobalInstance [static]

26.1.5.2 gdcmm::VRBINARY

## 26.2 gdcmm::network Namespace Reference

### Classes

- class [AAbortPDU](#)  
*[AAbortPDU Table](#) 9-26 A-ABORT PDU FIELDS.*
- class [AAssociateACPDU](#)  
*[AAssociateACPDU Table](#) 9-17 ASSOCIATE-AC PDU fields.*
- class [AAssociateRJPDU](#)  
*[AAssociateRJPDU Table](#) 9-21 ASSOCIATE-RJ PDU FIELDS.*
- class [AAssociateRQPDU](#)  
*[AAssociateRQPDU Table](#) 9-11 ASSOCIATE-RQ PDU fields.*
- class [AbstractSyntax](#)  
*[AbstractSyntax Table](#) 9-14 ABSTRACT SYNTAX SUB-ITEM FIELDS.*
- class [ApplicationContext](#)  
*[ApplicationContext Table](#) 9-12 APPLICATION CONTEXT ITEM FIELDS.*
- class [AReleaseRPPDU](#)  
*[AReleaseRPPDU Table](#) 9-25 A-RELEASE-RP PDU fields.*
- class [AReleaseRQPDU](#)  
*[AReleaseRQPDU Table](#) 9-24 A-RELEASE-RQ PDU FIELDS.*
- class [ARTIMTimer](#)  
*[ARTIMTimer](#) This file contains the code for the ARTIM timer.*
- class [AsynchronousOperationsWindowSub](#)  
*[AsynchronousOperationsWindowSub](#) PS 3.7 [Table D.3-7](#) ASYNCHRONOUS OPERATIONS WINDOW SUB-ITEM FIELDS (A-ASSOCIATE-RQ)↔*
- class [BaseCompositeMessage](#)  
*[BaseCompositeMessage](#) The Composite events described in section 3.7-2009 of the DICOM standard all use their own messages. These messages are constructed using Presentation Data Values, from section 3.8-2009 of the standard, and then fill in appropriate values in their datasets.*
- class [BasePDU](#)  
*[BasePDU](#) base class for PDUs.*
- class [CEchoRQ](#)  
*[CEchoRQ](#) this file defines the messages for the echo action.*
- class [CEchoRSP](#)  
*[CEchoRSP](#) this file defines the messages for the echo action.*
- class [CFind](#)
- class [CFindCancelRQ](#)  
*[CFindCancelRQ](#) this file defines the messages for the cfind action.*
- class [CFindRQ](#)  
*[CFindRQ](#) this file defines the messages for the cfind action.*

- class [CFindRSP](#)  
*CFindRSP* this file defines the messages for the cfind action.
- class [CMoveCancelRq](#)
- class [CMoveRQ](#)  
*CMoveRQ* this file defines the messages for the cmove action.
- class [CMoveRSP](#)  
*CMoveRSP* this file defines the messages for the cmove action.
- class [CompositeMessageFactory](#)  
*CompositeMessageFactory* This class constructs PDataPDUs, but that have been specifically constructed for the composite DICOM services (C-Echo, C-Find, C-Get, C-Move, and C-Store). It will also handle parsing the incoming data to determine which of the CompositePDUs the incoming data is, and so therefore allowing the scu to determine what to do with incoming data (if acting as a storescp server, for instance).
- class [CStoreRQ](#)  
*CStoreRQ* this file defines the messages for the cecho action.
- class [CStoreRSP](#)  
*CStoreRSP* this file defines the messages for the cecho action.
- class [DIMSE](#)  
*DIMSE* PS 3.7 - 2009 Annex E [Command](#) Dictionary (Normative) E.1 REGISTRY OF DICOM COMMAND ELEMENTS  
[Table E.1-1 COMMAND FIELDS \(PART 1\)](#)
- class [ImplementationClassUIDSub](#)  
*ImplementationClassUIDSub* PS 3.7 [Table D.3-1 IMPLEMENTATION CLASS UID SUB-ITEM FIELDS \(A-ASSOCIATE-RQ\)](#)
- class [ImplementationUIDSub](#)  
*ImplementationUIDSub* [Table D.3-2 IMPLEMENTATION UID SUB-ITEM FIELDS \(A-ASSOCIATE-AC\)](#)
- class [ImplementationVersionNameSub](#)  
*ImplementationVersionNameSub* [Table D.3-3 IMPLEMENTATION VERSION NAME SUB-ITEM FIELDS \(A-ASSOCIATE-RQ\)](#)
- class [MaximumLengthSub](#)  
*MaximumLengthSub* Annex D [Table D.1-1 MAXIMUM LENGTH SUB-ITEM FIELDS \(A-ASSOCIATE-RQ\)](#)
- class [PDataTFPDU](#)  
*PDataTFPDU* [Table 9-22 P-DATA-TF PDU FIELDS.](#)
- class [PDUFactory](#)  
*PDUFactory* basically, given an initial byte, construct the appropriate PDU. This way, the event loop doesn't have to know about all the different PDU types.
- class [PresentationContextAC](#)  
*PresentationContextAC* [Table 9-18 PRESENTATION CONTEXT ITEM FIELDS.](#)
- class [PresentationContextRQ](#)  
*PresentationContextRQ* [Table 9-13 PRESENTATION CONTEXT ITEM FIELDS.](#)
- class [PresentationDataValue](#)  
*PresentationDataValue* [Table 9-23 PRESENTATION-DATA-VALUE ITEM FIELDS.](#)
- class [RoleSelectionSub](#)  
*RoleSelectionSub* PS 3.7 [Table D.3-9 SCP/SCU ROLE SELECTION SUB-ITEM FIELDS \(A-ASSOCIATE-RQ\)](#)
- class [ServiceClassApplicationInformation](#)
- class [SOPClassExtendedNegociationSub](#)  
*SOPClassExtendedNegociationSub* PS 3.7 [Table D.3-11 SOP CLASS EXTENDED NEGOTIATION SUB-ITEM FIELDS \(A-ASSOCIATE-RQ and A-ASSOCIATE-AC\)](#)
- class [TableRow](#)
- class [TransferSyntaxSub](#)  
*TransferSyntaxSub* [Table 9-15 TRANSFER SYNTAX SUB-ITEM FIELDS.](#)



- struct [Transition](#)
- class [ULAction](#)

*[ULAction](#) A [ULConnection](#) in a given [ULState](#) can perform certain [ULActions](#). This base class provides the interface for running those [ULActions](#) on a given [ULConnection](#).*

- class [ULActionAA1](#)
- class [ULActionAA2](#)
- class [ULActionAA3](#)
- class [ULActionAA4](#)
- class [ULActionAA5](#)
- class [ULActionAA6](#)
- class [ULActionAA7](#)
- class [ULActionAA8](#)
- class [ULActionAE1](#)
- class [ULActionAE2](#)
- class [ULActionAE3](#)
- class [ULActionAE4](#)
- class [ULActionAE5](#)
- class [ULActionAE6](#)
- class [ULActionAE7](#)
- class [ULActionAE8](#)
- class [ULActionAR1](#)
- class [ULActionAR10](#)
- class [ULActionAR2](#)
- class [ULActionAR3](#)
- class [ULActionAR4](#)
- class [ULActionAR5](#)
- class [ULActionAR6](#)
- class [ULActionAR7](#)
- class [ULActionAR8](#)
- class [ULActionAR9](#)
- class [ULActionDT1](#)
- class [ULActionDT2](#)
- class [ULBasicCallback](#)

*[ULBasicCallback](#) This is the most basic of callbacks for how the [ULConnectionManager](#) handles incoming datasets. [DataSets](#) are just concatenated to the [mDataSets](#) vector, and the result can be pulled out of the vector by later code. Alternatives to this method include progress updates, saving to disk, etc. This class is NOT THREAD SAFE. Access the dataset vector after the entire set of datasets has been returned by the [ULConnectionManager](#).*

- class [ULConnection](#)

*[ULConnection](#) This is the class that contains the socket to another machine, and passes data through itself, as well as maintaining a sense of state.*

- class [ULConnectionCallback](#)
- class [ULConnectionInfo](#)

*[ULConnectionInfo](#) this class contains all the information about a particular connection as established by the user. That is, it's: User Information Calling AE Title Called AE Title IP address/computer name IP Port A connection must be established with this information, that's subsequently placed into various primitives for actual communication.*

- class [ULConnectionManager](#)

*[ULConnectionManager](#) The [ULConnectionManager](#) performs actions on the [ULConnection](#) given inputs from the user and from the state of what's going on around the connection (ie, timeouts of the ARTIM timer, responses from the peer across the connection, etc).*

- class [ULEvent](#)

*[ULEvent](#) base class for network events.*

- class [ULTransitionTable](#)

*ULTransitionTable* The transition table of all the ULEvents, new ULActions, and ULStates.

- class [ULWritingCallback](#)
- class [UserInformation](#)

*UserInformation Table 9-16 USER INFORMATION ITEM FIELDS.*

## Enumerations

- enum [EEventID](#) {  
[eAASSOCIATERequestLocalUser](#) = 0,  
[eTransportConnConfirmLocal](#),  
[eASSOCIATE\\_ACPDUreceived](#),  
[eASSOCIATE\\_RJPDUreceived](#),  
[eTransportConnIndicLocal](#),  
[eAASSOCIATE\\_RQPDUreceived](#),  
[eAASSOCIATEResponseAccept](#),  
[eAASSOCIATEResponseReject](#),  
[ePDATArequest](#),  
[ePDATATFPDU](#),  
[eARELEASERequest](#),  
[eARELEASE\\_RQPDUReceivedOpen](#),  
[eARELEASE\\_RPPDUReceived](#),  
[eARELEASEResponse](#),  
[eAABORTRequest](#),  
[eAABORTPDUReceivedOpen](#),  
[eTransportConnectionClosed](#),  
[eARTIMTimerExpired](#),  
[eUnrecognizedPDUReceived](#),  
[eEventDoesNotExist](#) }
- enum [EStateID](#) {  
[eStaDoesNotExist](#) = 0,  
[eSta1Idle](#) = 1,  
[eSta2Open](#) = 2,  
[eSta3WaitLocalAssoc](#) = 4,  
[eSta4LocalAssocDone](#) = 8,  
[eSta5WaitRemoteAssoc](#) = 16,  
[eSta6TransferReady](#) = 32,  
[eSta7WaitRelease](#) = 64,  
[eSta8WaitLocalRelease](#) = 128,  
[eSta9ReleaseCollisionRqLocal](#) = 256,  
[eSta10ReleaseCollisionAc](#) = 512,  
[eSta11ReleaseCollisionRq](#) = 1024,  
[eSta12ReleaseCollisionAcLocal](#) = 2048,  
[eSta13AwaitingClose](#) = 4096 }

## Functions

- int [GetStateIndex](#) ([EStateID](#) inState)

## Variables

- const int [cMaxEventID](#) = [eEventDoesNotExist](#)

- const int [cMaxStateID](#) = 13

## 26.2.1 Enumeration Type Documentation

### 26.2.1.1 enum gdcmm::network::EEventID

Enumerator

***eAASSOCIATERequestLocalUser***  
***eTransportConnConfirmLocal***  
***eASSOCIATE\_ACPDUreceived***  
***eASSOCIATE\_RJPDUreceived***  
***eTransportConnIndicLocal***  
***eAASSOCIATE\_RQPDUreceived***  
***eAASSOCIATEResponseAccept***  
***eAASSOCIATEResponseReject***  
***ePDATArequest***  
***ePDATATFPDU***  
***eARELEASERequest***  
***eARELEASE\_RQPDUReceivedOpen***  
***eARELEASE\_RPPDUReceived***  
***eARELEASEResponse***  
***eAABORTRequest***  
***eAABORTPDUReceivedOpen***  
***eTransportConnectionClosed***  
***eARTIMTimerExpired***  
***eUnrecognizedPDUReceived***  
***eEventDoesNotExist***

### 26.2.1.2 enum gdcmm::network::EStateID

Each network connection will be in a particular state at any given time. Those states have IDs as described in the standard ps3.8-2009, roughly 1-13. This enumeration lists those states. The actual ULState class will contain more information about transitions to other states.

name and date: 16 sept 2010 mmr

Enumerator

***eStaDoesNotExist***  
***eSta1Idle***  
***eSta2Open***  
***eSta3WaitLocalAssoc***  
***eSta4LocalAssocDone***  
***eSta5WaitRemoteAssoc***  
***eSta6TransferReady***

***eSta7WaitRelease***  
***eSta8WaitLocalRelease***  
***eSta9ReleaseCollisionRqLocal***  
***eSta10ReleaseCollisionAc***  
***eSta11ReleaseCollisionRq***  
***eSta12ReleaseCollisionAcLocal***  
***eSta13AwaitingClose***

## 26.2.2 Function Documentation

26.2.2.1 `int gdcn::network::GetStateIndex ( EStateID inState ) [inline]`

References `eSta10ReleaseCollisionAc`, `eSta11ReleaseCollisionRq`, `eSta12ReleaseCollisionAcLocal`, `eSta13AwaitingClose`, `eSta1Idle`, `eSta2Open`, `eSta3WaitLocalAssoc`, `eSta4LocalAssocDone`, `eSta5WaitRemoteAssoc`, `eSta6TransferReady`, `eSta7WaitRelease`, `eSta8WaitLocalRelease`, `eSta9ReleaseCollisionRqLocal`, and `eStaDoesNotExist`.

## 26.2.3 Variable Documentation

26.2.3.1 `const int gdcn::network::cMaxEventID = eEventDoesNotExist`

26.2.3.2 `const int gdcn::network::cMaxStateID = 13`

Referenced by `gdcn::network::TableRow::TableRow()`, and `gdcn::network::TableRow::~~TableRow()`.

## 26.3 gdcn::SegmentHelper Namespace Reference

### Classes

- struct [BasicCodedEntry](#)  
*This structure defines a basic coded entry with all of its attributes.*

## 26.4 gdcn::terminal Namespace Reference

Class for Terminal Allow one to print in color in a shell.

### Enumerations

- enum [Attribute](#) {  
`reset = 0,`  
`bright = 1,`  
`dim = 2,`  
`underline = 3,`  
`blink = 5,`  
`reverse = 7,`  
`hidden = 8` }

- enum `Color` {  
    `black` = 0,  
    `red`,  
    `green`,  
    `yellow`,  
    `blue`,  
    `magenta`,  
    `cyan`,  
    `white` }
- enum `Mode` {  
    `CONSOLE` = 0,  
    `VT100` }

## Functions

- `GDCM_EXPORT` std::string `setattribute` (`Attribute` att)
- `GDCM_EXPORT` std::string `setbgcolor` (`Color` c)
- `GDCM_EXPORT` std::string `setfgcolor` (`Color` c)
- `GDCM_EXPORT` void `setmode` (`Mode` m)

### 26.4.1 Detailed Description

Class for Terminal Allow one to print in color in a shell.

- support VT100 compatible shell
- win32 console

### 26.4.2 Enumeration Type Documentation

#### 26.4.2.1 enum `gdcmm::terminal::Attribute`

Enumerator

***reset***  
***bright***  
***dim***  
***underline***  
***blink***  
***reverse***  
***hidden***

#### 26.4.2.2 enum `gdcmm::terminal::Color`

Enumerator

***black***  
***red***  
***green***

*yellow*

*blue*

*magenta*

*cyan*

*white*

#### 26.4.2.3 enum `gdcmm::terminal::Mode`

Enumerator

***CONSOLE***

***VT100***

#### 26.4.3 Function Documentation

26.4.3.1 `GDCM_EXPORT std::string gdcmm::terminal::setattribute ( Attribute att )`

26.4.3.2 `GDCM_EXPORT std::string gdcmm::terminal::setbgcolor ( Color c )`

26.4.3.3 `GDCM_EXPORT std::string gdcmm::terminal::setfgcolor ( Color c )`

26.4.3.4 `GDCM_EXPORT void gdcmm::terminal::setmode ( Mode m )`

## Chapter 27


# Class Documentation

### 27.1 gdcm::network::AAabortPDU Class Reference


[AAabortPDU](#) [Table 9-26](#) A-ABORT PDU FIELDS.

```
#include <gdcmAAabortPDU.h>
```

Inheritance diagram for `gdcm::network::AAabortPDU`:

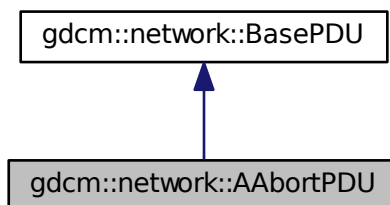

Collaboration diagram for `gdcm::network::AAabortPDU`:

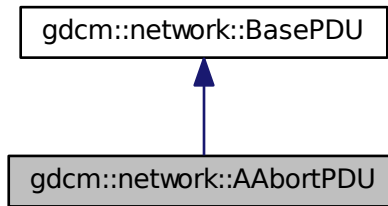

## Public Member Functions

- [AAabortPDU](#) ()
- `bool` [IsLastFragment](#) () const
- `void` [Print](#) (std::ostream &os) const
- `std::istream &` [Read](#) (std::istream &is)
- `void` [SetReason](#) (const uint8\_t r)
- `void` [SetSource](#) (const uint8\_t s)
- `size_t` [Size](#) () const
- `const std::ostream &` [Write](#) (std::ostream &os) const

### 27.1.1 Detailed Description

[AAabortPDU Table](#) 9-26 A-ABORT PDU FIELDS.

### 27.1.2 Constructor & Destructor Documentation

27.1.2.1 `gdcm::network::AAabortPDU::AAabortPDU ( )`

### 27.1.3 Member Function Documentation

27.1.3.1 `bool` `gdcm::network::AAabortPDU::IsLastFragment ( )` const `[inline]`,`[virtual]`

Implements [gdcm::network::BasePDU](#).

27.1.3.2 `void` `gdcm::network::AAabortPDU::Print ( std::ostream & os )` const `[virtual]`

Implements [gdcm::network::BasePDU](#).

27.1.3.3 `std::istream&` `gdcm::network::AAabortPDU::Read ( std::istream & is )` `[virtual]`

Implements [gdcm::network::BasePDU](#).



27.1.3.4 void gdcm::network::AAbortPDU::SetReason ( const uint8\_t r )

27.1.3.5 void gdcm::network::AAbortPDU::SetSource ( const uint8\_t s )

27.1.3.6 size\_t gdcm::network::AAbortPDU::Size ( ) const [virtual]

Implements [gdcm::network::BasePDU](#).

27.1.3.7 const std::ostream& gdcm::network::AAbortPDU::Write ( std::ostream & os ) const [virtual]

Implements [gdcm::network::BasePDU](#).

The documentation for this class was generated from the following file:


- [gdcmAAbortPDU.h](#)

## 27.2 gdcm::network::AAssociateACPDU Class Reference


[AAssociateACPDU](#) Table 9-17 ASSOCIATE-AC PDU fields.

```
#include <gdcmAAssociateACPDU.h>
```

Inheritance diagram for gdcm::network::AAssociateACPDU:

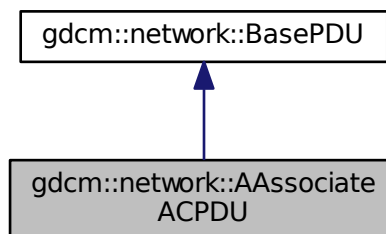

Collaboration diagram for `gdcmm::network::AAssociateACPDU`:

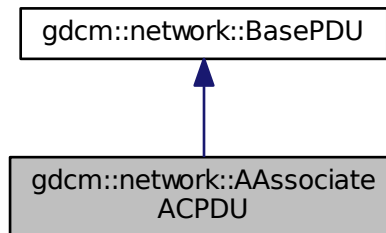

## Public Types

- typedef `std::vector< PresentationContextAC >::size_type` [SizeType](#)

## Public Member Functions

- [AAssociateACPDU](#) ()
- void [AddPresentationContextAC](#) ([PresentationContextAC](#) const &pcac)
- [SizeType](#) [GetNumberOfPresentationContextAC](#) () const
- const [PresentationContextAC](#) & [GetPresentationContextAC](#) ([SizeType](#) i)
- const [UserInformation](#) & [GetUserInformation](#) () const
- void [InitFromRQ](#) ([AAssociateRQPDU](#) const &rqpdu)
- bool [IsLastFragment](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- [SizeType](#) [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

## Protected Member Functions

- void [SetCalledAETitle](#) (const char calledaetitle[16])
- void [SetCallingAETitle](#) (const char callingaetitle[16])

## Friends

- class [AAssociateRQPDU](#)

## 27.2.1 Detailed Description

[AAssociateACPDU Table](#) 9-17 ASSOCIATE-AC PDU fields.

## 27.2.2 Member Typedef Documentation

27.2.2.1 `typedef std::vector<PresentationContextAC>::size_type gdcmm::network::AAAssociateACPDU::SizeType`

## 27.2.3 Constructor & Destructor Documentation

27.2.3.1 `gdcmm::network::AAAssociateACPDU::AAAssociateACPDU ( )`

## 27.2.4 Member Function Documentation

27.2.4.1 `void gdcmm::network::AAAssociateACPDU::AddPresentationContextAC ( PresentationContextAC const & pcac )`

27.2.4.2 `SizeType gdcmm::network::AAAssociateACPDU::GetNumberOfPresentationContextAC ( ) const [inline]`

27.2.4.3 `const PresentationContextAC& gdcmm::network::AAAssociateACPDU::GetPresentationContextAC ( SizeType i ) [inline]`

27.2.4.4 `const UserInformation& gdcmm::network::AAAssociateACPDU::GetUserInformation ( ) const [inline]`

27.2.4.5 `void gdcmm::network::AAAssociateACPDU::InitFromRQ ( AAAssociateRQPDU const & rqpdu )`

27.2.4.6 `bool gdcmm::network::AAAssociateACPDU::IsLastFragment ( ) const [inline],[virtual]`

Implements [gdcmm::network::BasePDU](#).

27.2.4.7 `void gdcmm::network::AAAssociateACPDU::Print ( std::ostream & os ) const [virtual]`

Implements [gdcmm::network::BasePDU](#).

27.2.4.8 `std::istream& gdcmm::network::AAAssociateACPDU::Read ( std::istream & is ) [virtual]`

Implements [gdcmm::network::BasePDU](#).

27.2.4.9 `void gdcmm::network::AAAssociateACPDU::SetCalledAETitle ( const char calledaetitle[16] ) [protected]`

27.2.4.10 `void gdcmm::network::AAAssociateACPDU::SetCallingAETitle ( const char callingaetitle[16] ) [protected]`

27.2.4.11 `SizeType gdcmm::network::AAAssociateACPDU::Size ( ) const [virtual]`

Implements [gdcmm::network::BasePDU](#).

27.2.4.12 `const std::ostream& gdcmm::network::AAAssociateACPDU::Write ( std::ostream & os ) const [virtual]`

Implements [gdcmm::network::BasePDU](#).

## 27.2.5 Friends And Related Function Documentation

### 27.2.5.1 friend class **AAssociateRQPDU** [friend]

The documentation for this class was generated from the following file:


- [gdcmAAssociateACPDU.h](#)

## 27.3 gdcmm::network::AAssociateRJPDU Class Reference


[AAssociateRJPDU](#) Table 9-21 ASSOCIATE-RJ PDU FIELDS.

```
#include <gdcmAAssociateRJPDU.h>
```

Inheritance diagram for gdcmm::network::AAssociateRJPDU:

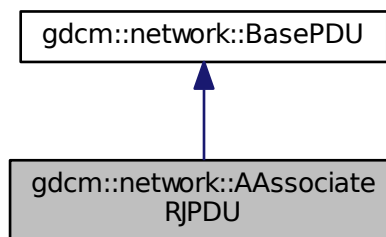

Collaboration diagram for gdcmm::network::AAssociateRJPDU:

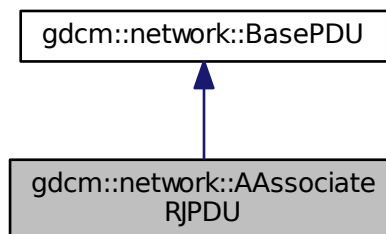

### Public Member Functions

- [AAssociateRJPDU](#) ()
- bool [IsLastFragment](#) () const

- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.3.1 Detailed Description

[AAssociateRJPDUTable](#) 9-21 ASSOCIATE-RJ PDU FIELDS.

### 27.3.2 Constructor & Destructor Documentation

27.3.2.1 `gdcm::network::AAssociateRJPDUTable::AAssociateRJPDUTable ( )`

### 27.3.3 Member Function Documentation

27.3.3.1 `bool gdcm::network::AAssociateRJPDUTable::IsLastFragment ( ) const` `[inline]`, `[virtual]`

Implements [gdcm::network::BasePDU](#).

27.3.3.2 `void gdcm::network::AAssociateRJPDUTable::Print ( std::ostream & os ) const` `[virtual]`

Implements [gdcm::network::BasePDU](#).

27.3.3.3 `std::istream& gdcm::network::AAssociateRJPDUTable::Read ( std::istream & is )` `[virtual]`

Implements [gdcm::network::BasePDU](#).

27.3.3.4 `size_t gdcm::network::AAssociateRJPDUTable::Size ( ) const` `[virtual]`

Implements [gdcm::network::BasePDU](#).

27.3.3.5 `const std::ostream& gdcm::network::AAssociateRJPDUTable::Write ( std::ostream & os ) const` `[virtual]`

Implements [gdcm::network::BasePDU](#).

The documentation for this class was generated from the following file:


- [gdcmAAssociateRJPDUTable.h](#)

## 27.4 gdcm::network::AAssociateRQPDU Class Reference


[AAssociateRQPDU Table](#) 9-11 ASSOCIATE-RQ PDU fields.

```
#include <gdcmAAssociateRQPDU.h>
```

Inheritance diagram for `gdc::network::AAssociateRQPDU`:

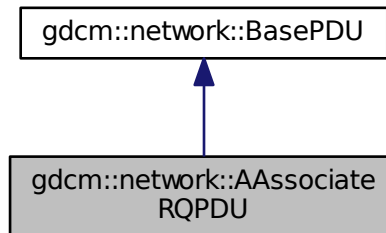

Collaboration diagram for `gdc::network::AAssociateRQPDU`:

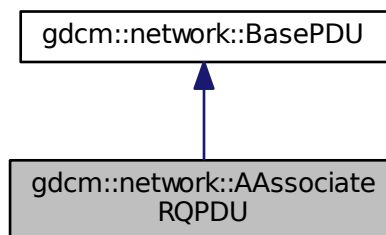

## Public Types

- `typedef std::vector< PresentationContextRQ > PresentationContextArrayType`
- `typedef std::vector< PresentationContextRQ >::size_type SizeType`

## Public Member Functions

- [AAssociateRQPDU](#) ()
- [AAssociateRQPDU](#) (const [AAssociateRQPDU](#) &pdu)
- void [AddPresentationContext](#) ([PresentationContextRQ](#) const &pc)
- std::string [GetCalledAETitle](#) () const
- std::string [GetCallingAETitle](#) () const
- [SizeType](#) [GetNumberOfPresentationContext](#) () const
- [PresentationContextRQ](#) const & [GetPresentationContext](#) ([SizeType](#) i) const
- const [PresentationContextRQ](#) \* [GetPresentationContextByAbstractSyntax](#) ([AbstractSyntax](#) const &as) const
- const [PresentationContextRQ](#) \* [GetPresentationContextByID](#) (uint8\_t i) const

- [PresentationContextArrayType](#) const & [GetPresentationContexts](#) ()
- const [UserInformation](#) & [GetUserInformation](#) () const
- bool [IsLastFragment](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetCalledAETitle](#) (const char calledaetitle[16])  
*Set the Called AE Title.*
- void [SetCallingAETitle](#) (const char callingaetitle[16])  
*Set the Calling AE Title.*
- void [SetUserInformation](#) ([UserInformation](#) const &ui)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### Static Public Member Functions

- static bool [IsAETitleValid](#) (const char title[16])  
*Check whether or not the.*

### Protected Member Functions

- std::string [GetReserved43\\_74](#) () const

### Friends

- class [AAssociateACPDU](#)

## 27.4.1 Detailed Description

[AAssociateRQPDU](#) Table 9-11 ASSOCIATE-RQ PDU fields.

## 27.4.2 Member Typedef Documentation

27.4.2.1 typedef std::vector<[PresentationContextRQ](#)> gdcm::network::AAssociateRQPDU::PresentationContextArrayType

27.4.2.2 typedef std::vector<[PresentationContextRQ](#)>::size\_type gdcm::network::AAssociateRQPDU::SizeType

## 27.4.3 Constructor & Destructor Documentation

27.4.3.1 gdcm::network::AAssociateRQPDU::AAssociateRQPDU ( )

27.4.3.2 gdcm::network::AAssociateRQPDU::AAssociateRQPDU ( const AAssociateRQPDU & pdu ) [inline]

## 27.4.4 Member Function Documentation

27.4.4.1 void gdcm::network::AAssociateRQPDU::AddPresentationContext ( [PresentationContextRQ](#) const & pc )

- 27.4.4.2 `std::string gdcn::network::AAAssociateRQPDU::GetCalledAETitle ( ) const` `[inline]`
- 27.4.4.3 `std::string gdcn::network::AAAssociateRQPDU::GetCallingAETitle ( ) const` `[inline]`
- 27.4.4.4 `SizeType gdcn::network::AAAssociateRQPDU::GetNumberOfPresentationContext ( ) const` `[inline]`
- 27.4.4.5 `PresentationContextRQ const& gdcn::network::AAAssociateRQPDU::GetPresentationContext ( SizeType i ) const` `[inline]`
- 27.4.4.6 `const PresentationContextRQ* gdcn::network::AAAssociateRQPDU::GetPresentationContextByAbstractSyntax ( AbstractSyntax const & as ) const`
- 27.4.4.7 `const PresentationContextRQ* gdcn::network::AAAssociateRQPDU::GetPresentationContextByID ( uint8_t i ) const`
- 27.4.4.8 `PresentationContextArrayType const& gdcn::network::AAAssociateRQPDU::GetPresentationContexts ( )` `[inline]`
- 27.4.4.9 `std::string gdcn::network::AAAssociateRQPDU::GetReserved43_74 ( ) const` `[protected]`
- 27.4.4.10 `const UserInformation& gdcn::network::AAAssociateRQPDU::GetUserInformation ( ) const` `[inline]`
- 27.4.4.11 `static bool gdcn::network::AAAssociateRQPDU::IsAETitleValid ( const char title[16] )` `[static]`

Check whether or not the.

Parameters

<i>title</i>	is a valid AE title
--------------	---------------------

- 27.4.4.12 `bool gdcn::network::AAAssociateRQPDU::IsLastFragment ( ) const` `[inline]`, `[virtual]`

Implements [gdcn::network::BasePDU](#).

- 27.4.4.13 `void gdcn::network::AAAssociateRQPDU::Print ( std::ostream & os ) const` `[virtual]`

This function will initialize an [AAAssociateACPDU](#) from the fields in the [AAAssociateRQPDU](#) structure

Implements [gdcn::network::BasePDU](#).

- 27.4.4.14 `std::istream& gdcn::network::AAAssociateRQPDU::Read ( std::istream & is )` `[virtual]`

Implements [gdcn::network::BasePDU](#).

- 27.4.4.15 `void gdcn::network::AAAssociateRQPDU::SetCalledAETitle ( const char calledaetitle[16] )`

Set the Called AE Title.

- 27.4.4.16 `void gdcn::network::AAAssociateRQPDU::SetCallingAETitle ( const char callingaetitle[16] )`

Set the Calling AE Title.



27.4.4.17 `void gdcm::network::AAssociateRQPDU::SetUserInfo ( UserInfo const & ui )`

27.4.4.18 `size_t gdcm::network::AAssociateRQPDU::Size ( ) const` [virtual]

Implements [gdcm::network::BasePDU](#).

27.4.4.19 `const std::ostream& gdcm::network::AAssociateRQPDU::Write ( std::ostream & os ) const` [virtual]

Implements [gdcm::network::BasePDU](#).

## 27.4.5 Friends And Related Function Documentation

27.4.5.1 `friend class AAssociateACPDU` [friend]


The documentation for this class was generated from the following file:

- [gdcmAAssociateRQPDU.h](#)


## 27.5 gdcm::AbortEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::AbortEvent`:

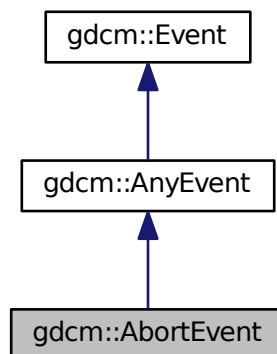

Collaboration diagram for `gdcm::AbortEvent`:

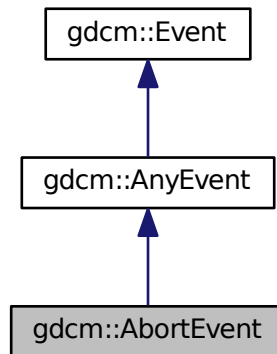

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.6 `gdcm::network::AbstractSyntax` Class Reference

[AbstractSyntax](#) Table 9-14 ABSTRACT SYNTAX SUB-ITEM FIELDS.

```
#include <gdcmAbstractSyntax.h>
```

### Public Member Functions

- [AbstractSyntax](#) ()
- [DataElement GetAsDataElement](#) () const
- const char \* [GetName](#) () const
- bool [operator==](#) (const [AbstractSyntax](#) &as) const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetName](#) (const char \*name)
- void [SetNameFromUID](#) (UIDs::TSName tsname)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

#### 27.6.1 Detailed Description

[AbstractSyntax](#) Table 9-14 ABSTRACT SYNTAX SUB-ITEM FIELDS.

## 27.6.2 Constructor & Destructor Documentation

27.6.2.1 `gdcm::network::AbstractSyntax::AbstractSyntax ( )`

## 27.6.3 Member Function Documentation

27.6.3.1 `DataElement gdcm::network::AbstractSyntax::GetAsDataElement ( ) const`

27.6.3.2 `const char* gdcm::network::AbstractSyntax::GetName ( ) const` `[inline]`

27.6.3.3 `bool gdcm::network::AbstractSyntax::operator== ( const AbstractSyntax & as ) const` `[inline]`

27.6.3.4 `void gdcm::network::AbstractSyntax::Print ( std::ostream & os ) const`

27.6.3.5 `std::istream& gdcm::network::AbstractSyntax::Read ( std::istream & is )`

27.6.3.6 `void gdcm::network::AbstractSyntax::SetName ( const char * name )` `[inline]`

27.6.3.7 `void gdcm::network::AbstractSyntax::SetNameFromUID ( UIDs::TSName tsname )`

27.6.3.8 `size_t gdcm::network::AbstractSyntax::Size ( ) const`

27.6.3.9 `const std::ostream& gdcm::network::AbstractSyntax::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:


- [gdcmAbstractSyntax.h](#)

## 27.7 gdcm::AnonymizeEvent Class Reference


[AnonymizeEvent](#) Special type of event triggered during the Anonymization process.

```
#include <gdcmAnonymizeEvent.h>
```

Inheritance diagram for `gdcm::AnonymizeEvent`:

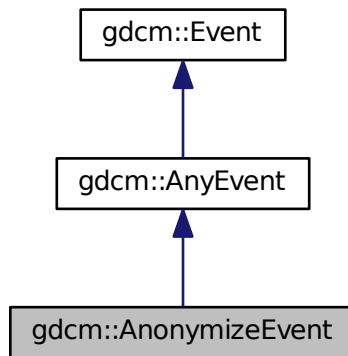

Collaboration diagram for `gdcm::AnonymizeEvent`:

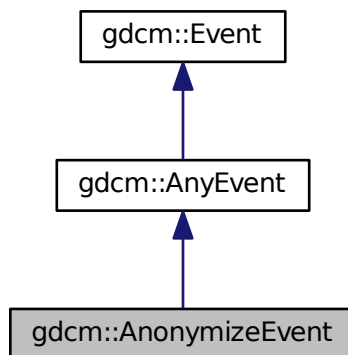

## Public Types

- typedef [AnonymizeEvent](#) `Self`
- typedef [AnyEvent](#) `Superclass`

## Public Member Functions

- [AnonymizeEvent](#) (`Tag` const &tag=0)
- [AnonymizeEvent](#) (const [Self](#) &s)

- virtual [~AnonymizeEvent](#) ()
- virtual bool [CheckEvent](#) (const [::gdcm::Event](#) \*e) const
- virtual const char \* [GetEventName](#) () const
- [Tag](#) const & [GetTag](#) () const
- virtual [::gdcm::Event](#) \* [MakeObject](#) () const
- void [SetTag](#) (const [Tag](#) &t)

### 27.7.1 Detailed Description

[AnonymizeEvent](#) Special type of event triggered during the Anonymization process.

See also

[Anonymizer](#)

### 27.7.2 Member Typedef Documentation

27.7.2.1 `typedef AnonymizeEvent gdcm::AnonymizeEvent::Self`

27.7.2.2 `typedef AnyEvent gdcm::AnonymizeEvent::Superclass`

### 27.7.3 Constructor & Destructor Documentation

27.7.3.1 `gdcm::AnonymizeEvent::AnonymizeEvent ( Tag const & tag = 0 ) [inline]`

27.7.3.2 `virtual gdcm::AnonymizeEvent::~~AnonymizeEvent ( ) [inline],[virtual]`

27.7.3.3 `gdcm::AnonymizeEvent::AnonymizeEvent ( const Self & s ) [inline]`

### 27.7.4 Member Function Documentation

27.7.4.1 `virtual bool gdcm::AnonymizeEvent::CheckEvent ( const ::gdcm::Event * e ) const [inline],[virtual]`

27.7.4.2 `virtual const char* gdcm::AnonymizeEvent::GetEventName ( ) const [inline],[virtual]`

Return the StringName associated with the event.

Implements [gdcm::Event](#).

27.7.4.3 `Tag const& gdcm::AnonymizeEvent::GetTag ( ) const [inline]`

27.7.4.4 `virtual ::gdcm::Event* gdcm::AnonymizeEvent::MakeObject ( ) const [inline],[virtual]`

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implements [gdcm::Event](#).

27.7.4.5 `void gdcm::AnonymizeEvent::SetTag ( const Tag & t ) [inline]`

The documentation for this class was generated from the following file:


- [gdcmAnonymizeEvent.h](#)

## 27.8 gdcmm::Anonymizer Class Reference


[Anonymizer](#) This class is a multi purpose anonymizer. It can work in 2 mode:

```
#include <gdcmmAnonymizer.h>
```

Inheritance diagram for gdcmm::Anonymizer:

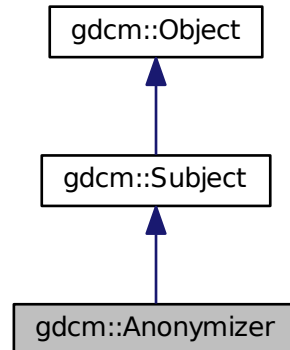

Collaboration diagram for gdcmm::Anonymizer:

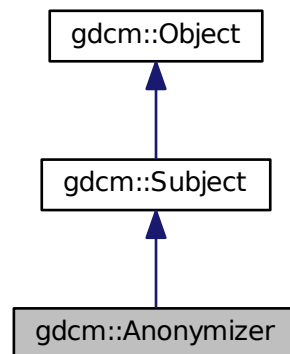

### Public Member Functions

- [Anonymizer](#) ()
- [~Anonymizer](#) ()

- bool [BasicApplicationLevelConfidentialityProfile](#) (bool deidentify=true)
- bool [Empty](#) (Tag const &t)
- const [CryptographicMessageSyntax](#) \* [GetCryptographicMessageSyntax](#) () const
- [File](#) & [GetFile](#) ()
- bool [Remove](#) (Tag const &t)
- bool [RemoveGroupLength](#) ()  
*Main function that loop over all elements and remove group length.*
- bool [RemovePrivateTags](#) ()  
*Main function that loop over all elements and remove private tags.*
- bool [RemoveRetired](#) ()  
*Main function that loop over all elements and remove retired element.*
- bool [Replace](#) (Tag const &t, const char \*value)
- bool [Replace](#) (Tag const &t, const char \*value, [VL](#) const &vl)
- void [SetCryptographicMessageSyntax](#) ([CryptographicMessageSyntax](#) \*cms)  
*Set/Get CMS key that will be used to encrypt the dataset within BasicApplicationLevelConfidentialityProfile.*
- void [SetFile](#) (const [File](#) &f)  
*Set/Get File.*

### Static Public Member Functions

- static void [ClearInternalUIDs](#) ()
- static std::vector< [Tag](#) > [GetBasicApplicationLevelConfidentialityProfileAttributes](#) ()  
*Return the list of Tag that will be considered when anonymizing a DICOM file.*
- static [SmartPointer](#)< [Anonymizer](#) > [New](#) ()  
*for wrapped language: instantiate a reference counted object*

### Protected Member Functions

- bool [BALCPPProtect](#) ([DataSet](#) &ds, [Tag](#) const &tag, const [IOD](#) &iod)
- bool [CanEmptyTag](#) (Tag const &tag, const [IOD](#) &iod) const
- void [RecurseDataSet](#) ([DataSet](#) &ds)

#### 27.8.1 Detailed Description

[Anonymizer](#) This class is a multi purpose anonymizer. It can work in 2 mode:

- Full (irreversible) anonymizer (aka dumb mode)
- reversible de-identifier/re-identifier (aka smart mode). This implements the Basic Application Level Confidentiality Profile, DICOM PS 3.15-2009

1. dumb mode This is a dumb anonymizer implementation. All it allows user is simple operation such as:

[Tag](#) based functions:

- complete removal of DICOM attribute (Remove)
- make a tag empty, ie make it's length 0 (Empty)

- replace with another string-based value (Replace)

[DataSet](#) based functions:

- Remove all group length attribute from a DICOM dataset (Group Length element are deprecated, DICOM 2008)
- Remove all private attributes
- Remove all retired attributes

All function calls actually execute the user specified request. Previous implementation were calling a general Anonymize function but traversing a `std::set` is  $O(n)$  operation, while a simple user specified request is  $O(\log(n))$  operation. So 'm' user interaction is  $O(m*\log(n))$  which is  $< O(n)$  complexity.

1. smart mode this mode implements the Basic Application Level Confidentiality Profile (DICOM PS 3.15-2008) In this case, it is extremely important to use the same [Anonymizer](#) class when anonymizing a [FileSet](#). Once the [Anonymizer](#) is destroyed its memory of known (already processed) [UIDs](#) will be lost. which will make the anonymizer behaves incorrectly for attributes such as [Series](#) [UID](#) [Study](#) [UID](#) where user want some consistency. When attribute is [Type 1](#) / [Type 1C](#), a dummy generator will take in the existing value and produce a dummy value (a sha1 representation). sha1 algorithm is considered to be cryptographically strong (compared to md5sum) so that we meet the following two conditions:

- Produce the same dummy value for the same input value
- do not provide an easy way to retrieve the original value from the sha1 generated value

This class implement the Subject/Observer pattern trigger the following event:

- [AnonymizeEvent](#)
- [IterationEvent](#)
- [StartEvent](#)
- [EndEvent](#)

See also

[CryptographicMessageSyntax](#)

Examples:

[ClinicalTrialAnnotate.cxx](#), [CreateJPIPDataSet.cxx](#), and [EncapsulateFileInRawData.cxx](#).

## 27.8.2 Constructor & Destructor Documentation

27.8.2.1 `gdcm::Anonymizer::Anonymizer ( )` [`inline`]

27.8.2.2 `gdcm::Anonymizer::~~Anonymizer ( )`

## 27.8.3 Member Function Documentation

27.8.3.1 `bool gdcm::Anonymizer::BALCPPProtect ( DataSet & ds, Tag const & tag, const IOD & iod )` [`protected`]

27.8.3.2 `bool gdcm::Anonymizer::BasicApplicationLevelConfidentialityProfile ( bool deidentify = true )`

PS 3.15 / E.1.1 De-Identifier An Application may claim conformance to the Basic Application Level Confidentiality Profile as a deidentifier if it protects all Attributes that might be used by unauthorized entities to identify the patient. NOT THREAD SAFE



27.8.3.3 `bool gdcm::Anonymizer::CanEmptyTag ( Tag const & tag, const IOD & iod ) const` `[protected]`

27.8.3.4 `static void gdcm::Anonymizer::ClearInternalUIDs ( )` `[static]`

Clear the internal mapping of real [UIDs](#) to generated [UIDs](#)

Warning

the mapping is definitely lost

27.8.3.5 `bool gdcm::Anonymizer::Empty ( Tag const & t )`

Make [Tag](#) t empty (if not found tag will be created) Warning: does not handle SQ element

Examples:

[CreateJPIPDataSet.cxx](#).

27.8.3.6 `static std::vector<Tag> gdcm::Anonymizer::GetBasicApplicationLevelConfidentialityProfileAttributes ( )` `[static]`

Return the list of [Tag](#) that will be considered when anonymizing a DICOM file.

Examples:

[GenFakeIdentifyFile.cxx](#), and [TraverseModules.cxx](#).

27.8.3.7 `const CryptographicMessageSyntax* gdcm::Anonymizer::GetCryptographicMessageSyntax ( ) const`

27.8.3.8 `File& gdcm::Anonymizer::GetFile ( )` `[inline]`

27.8.3.9 `static SmartPointer<Anonymizer> gdcm::Anonymizer::New ( )` `[inline], [static]`

for wrapped language: instantiate a reference counted object

27.8.3.10 `void gdcm::Anonymizer::RecurseDataSet ( DataSet & ds )` `[protected]`

27.8.3.11 `bool gdcm::Anonymizer::Remove ( Tag const & t )`

remove a tag (even a SQ can be removed) Return code is false when tag t cannot be found

27.8.3.12 `bool gdcm::Anonymizer::RemoveGroupLength ( )`

Main function that loop over all elements and remove group length.

Examples:

[ClinicalTrialAnnotate.cxx](#).

### 27.8.3.13 `bool gdcm::Anonymizer::RemovePrivateTags ( )`

Main function that loop over all elements and remove private tags.

Examples:

[ClinicalTrialAnnotate.cxx](#).

### 27.8.3.14 `bool gdcm::Anonymizer::RemoveRetired ( )`

Main function that loop over all elements and remove retired element.

### 27.8.3.15 `bool gdcm::Anonymizer::Replace ( Tag const & t, const char * value )`

Replace tag with another value, if tag is not found it will be created: WARNING: this function can only execute if tag is a VRASCI

Examples:

[ClinicalTrialAnnotate.cxx](#), [CreateJPIPDataSet.cxx](#), and [EncapsulateFileInRawData.cxx](#).

### 27.8.3.16 `bool gdcm::Anonymizer::Replace ( Tag const & t, const char * value, VL const & vl )`

when the value contains \0, it is a good idea to specify the length. This function is required when dealing with VRBINARY tag

### 27.8.3.17 `void gdcm::Anonymizer::SetCryptographicMessageSyntax ( CryptographicMessageSyntax * cms )`

Set/Get CMS key that will be used to encrypt the dataset within BasicApplicationLevelConfidentialityProfile.

### 27.8.3.18 `void gdcm::Anonymizer::SetFile ( const File & f ) [inline]`

Set/Get [File](#).

Examples:

[ClinicalTrialAnnotate.cxx](#), [CreateJPIPDataSet.cxx](#), and [EncapsulateFileInRawData.cxx](#).


The documentation for this class was generated from the following file:

- [gdcmAnonymizer.h](#)


## 27.9 `gdcm::AnyEvent` Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for gdcM::AnyEvent:

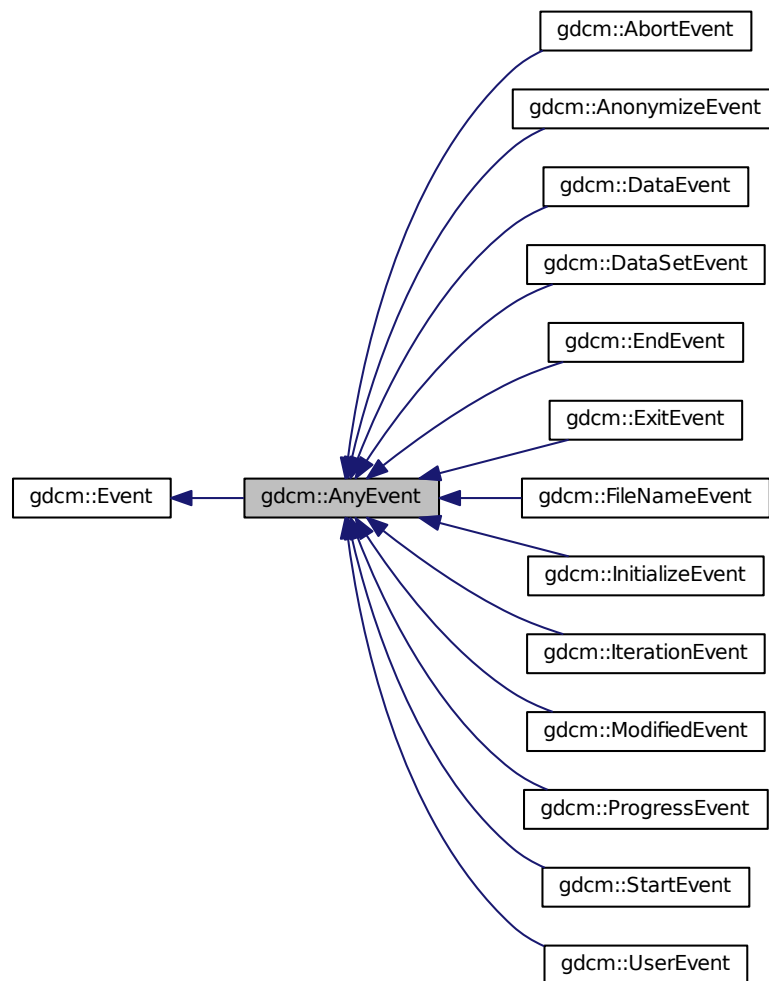

Collaboration diagram for `gdcm::AnyEvent`:

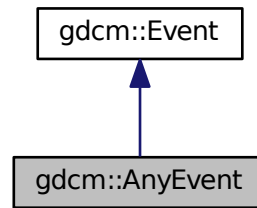

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.10 `gdcm::network::ApplicationContext` Class Reference

[ApplicationContext](#) Table 9-12 APPLICATION CONTEXT ITEM FIELDS.

```
#include <gdcmApplicationContext.h>
```

### Public Member Functions

- [ApplicationContext](#) ()
- `const char *` [GetName](#) () `const`
- `void` [Print](#) (`std::ostream &os`) `const`
- `std::istream &` [Read](#) (`std::istream &is`)
- `void` [SetName](#) (`const char *name`)
- `size_t` [Size](#) () `const`
- `const std::ostream &` [Write](#) (`std::ostream &os`) `const`

#### 27.10.1 Detailed Description

[ApplicationContext](#) Table 9-12 APPLICATION CONTEXT ITEM FIELDS.

**Todo** Looks like Application Context can only be 64 bytes at max (see Figure 9-1 / PS 3.8 - 2009 )

#### 27.10.2 Constructor & Destructor Documentation

27.10.2.1 `gdcm::network::ApplicationContext::ApplicationContext ( )`

### 27.10.3 Member Function Documentation

27.10.3.1 `const char* gdcm::network::ApplicationContext::GetName ( ) const` `[inline]`

27.10.3.2 `void gdcm::network::ApplicationContext::Print ( std::ostream & os ) const`

27.10.3.3 `std::istream& gdcm::network::ApplicationContext::Read ( std::istream & is )`

27.10.3.4 `void gdcm::network::ApplicationContext::SetName ( const char * name )` `[inline]`

27.10.3.5 `size_t gdcm::network::ApplicationContext::Size ( ) const`

27.10.3.6 `const std::ostream& gdcm::network::ApplicationContext::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:


- [gdcmApplicationContext.h](#)

## 27.11 gdcm::ApplicationEntity Class Reference

[ApplicationEntity](#).

```
#include <gdcmApplicationEntity.h>
```

Collaboration diagram for gdcm::ApplicationEntity:

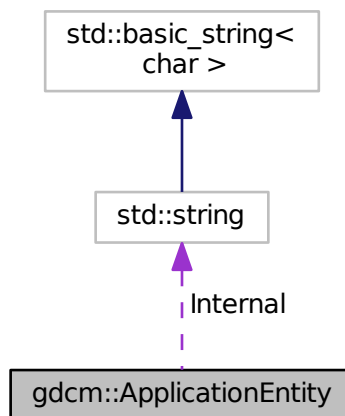

### Public Member Functions

- `bool IsValid ( ) const`
- `void Print (std::ostream & os) const`

- void [SetBlob](#) (const std::vector< char > &v)
- void [Squeeze](#) ()

### Public Attributes

- std::string [Internal](#)

### Static Public Attributes

- static const unsigned int [MaxLength](#) = 16
- static const unsigned int [MaxNumberOfComponents](#) = 1
- static const char [Padding](#) = ''
- static const char [Separator](#) = ''

## 27.11.1 Detailed Description

[ApplicationEntity](#).

- AE Application Entity
- A string of characters that identifies an Application Entity with leading and trailing spaces (20H) being non-significant. A value consisting solely of spaces shall not be used.
- Default Character Repertoire excluding character code 5CH (the BACKSLASH \ in ISO-IR 6), and control characters LF, FF, CR and ESC.
- 16 bytes maximum

## 27.11.2 Member Function Documentation

27.11.2.1 bool [gdcm::ApplicationEntity::IsValid](#) ( ) const [inline]

27.11.2.2 void [gdcm::ApplicationEntity::Print](#) ( std::ostream & os ) const [inline]

27.11.2.3 void [gdcm::ApplicationEntity::SetBlob](#) ( const std::vector< char > & v ) [inline]

27.11.2.4 void [gdcm::ApplicationEntity::Squeeze](#) ( ) [inline]

## 27.11.3 Member Data Documentation

27.11.3.1 std::string [gdcm::ApplicationEntity::Internal](#)

27.11.3.2 const unsigned int [gdcm::ApplicationEntity::MaxLength](#) = 16 [static]

27.11.3.3 const unsigned int [gdcm::ApplicationEntity::MaxNumberOfComponents](#) = 1 [static]

27.11.3.4 const char [gdcm::ApplicationEntity::Padding](#) = '' [static]

27.11.3.5 const char [gdcm::ApplicationEntity::Separator](#) = '' [static]

The documentation for this class was generated from the following file:


- [gdcmmApplicationEntity.h](#)

## 27.12 gdcmm::network::AReleaseRPPDU Class Reference


[AReleaseRPPDU](#) Table 9-25 A-RELEASE-RP PDU fields.

```
#include <gdcmmAReleaseRPPDU.h>
```

Inheritance diagram for gdcmm::network::AReleaseRPPDU:

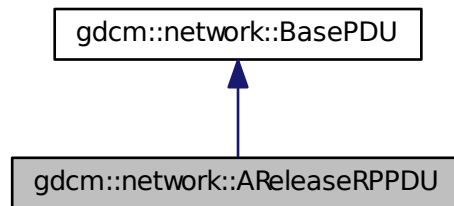

Collaboration diagram for gdcmm::network::AReleaseRPPDU:

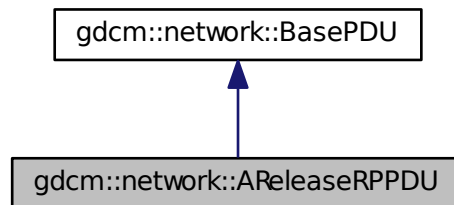

### Public Member Functions

- [AReleaseRPPDU](#) ()
- bool [IsLastFragment](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.12.1 Detailed Description

[AReleaseRPPDU Table](#) 9-25 A-RELEASE-RP PDU fields.

### 27.12.2 Constructor & Destructor Documentation

27.12.2.1 `gdcn::network::AReleaseRPPDU::AReleaseRPPDU ( )`

### 27.12.3 Member Function Documentation

27.12.3.1 `bool gdcn::network::AReleaseRPPDU::IsLastFragment ( ) const` `[inline],[virtual]`

Implements [gdcn::network::BasePDU](#).

27.12.3.2 `void gdcn::network::AReleaseRPPDU::Print ( std::ostream & os ) const` `[virtual]`

Implements [gdcn::network::BasePDU](#).

27.12.3.3 `std::istream& gdcn::network::AReleaseRPPDU::Read ( std::istream & is )` `[virtual]`

Implements [gdcn::network::BasePDU](#).

27.12.3.4 `size_t gdcn::network::AReleaseRPPDU::Size ( ) const` `[virtual]`

Implements [gdcn::network::BasePDU](#).

27.12.3.5 `const std::ostream& gdcn::network::AReleaseRPPDU::Write ( std::ostream & os ) const` `[virtual]`

Implements [gdcn::network::BasePDU](#).

The documentation for this class was generated from the following file:

- [gdcnAReleaseRPPDU.h](#)


## 27.13 gdcn::network::AReleaseRQPDU Class Reference

[AReleaseRQPDU Table](#) 9-24 A-RELEASE-RQ PDU FIELDS.


```
#include <gdcnAReleaseRQPDU.h>
```



Inheritance diagram for gdcmm::network::AReleaseRQPDU:

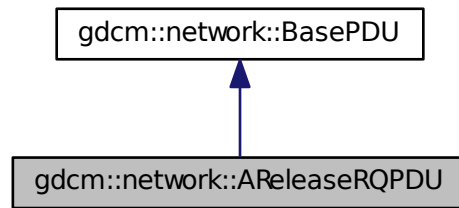

Collaboration diagram for gdcmm::network::AReleaseRQPDU:

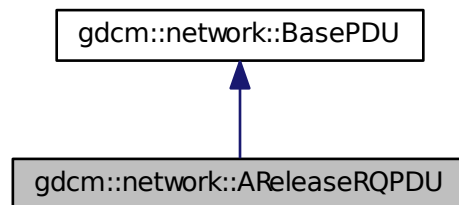

## Public Member Functions

- [AReleaseRQPDU](#) ()
- bool [IsLastFragment](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.13.1 Detailed Description

[AReleaseRQPDU](#) Table 9-24 A-RELEASE-RQ PDU FIELDS.

### 27.13.2 Constructor & Destructor Documentation

27.13.2.1 gdcmm::network::AReleaseRQPDU::AReleaseRQPDU ( )

### 27.13.3 Member Function Documentation

27.13.3.1 `bool gdcn::network::AReleaseRQPDU::IsLastFragment ( ) const [inline],[virtual]`

Implements [gdcn::network::BasePDU](#).

27.13.3.2 `void gdcn::network::AReleaseRQPDU::Print ( std::ostream & os ) const [virtual]`

Implements [gdcn::network::BasePDU](#).

27.13.3.3 `std::istream& gdcn::network::AReleaseRQPDU::Read ( std::istream & is ) [virtual]`

Implements [gdcn::network::BasePDU](#).

27.13.3.4 `size_t gdcn::network::AReleaseRQPDU::Size ( ) const [virtual]`

Implements [gdcn::network::BasePDU](#).

27.13.3.5 `const std::ostream& gdcn::network::AReleaseRQPDU::Write ( std::ostream & os ) const [virtual]`

Implements [gdcn::network::BasePDU](#).

The documentation for this class was generated from the following file:

- [gdcnAReleaseRQPDU.h](#)

## 27.14 gdcn::network::ARTIMTimer Class Reference

[ARTIMTimer](#) This file contains the code for the ARTIM timer.

```
#include <gdcnARTIMTimer.h>
```

### Public Member Functions

- [ARTIMTimer](#) ()
- double [GetElapsedTime](#) () const
- bool [GetHasExpired](#) () const
- double [GetTimeout](#) () const
- void [SetTimeout](#) (double inTimeout)
- void [Start](#) ()
- void [Stop](#) ()

### 27.14.1 Detailed Description

[ARTIMTimer](#) This file contains the code for the ARTIM timer.

Basically, the ARTIM timer will just get the wall time when it's started, and then can be queried for the current time, and then can be stopped (ie, the start time reset).

Because we're trying to do this without threading, we should be able to 'start' the ARTIM timer by this mechanism, and then when waiting for a particular response, tight loop that with sleep calls and determinations of when the ARTIM timer has reached its peak. As such, this isn't a strict 'timer' in the traditional sense of the word, but more of a time keeper.

There can be only one ARTIM timer per connection.

## 27.14.2 Constructor & Destructor Documentation

27.14.2.1 `gdcm::network::ARTIMTimer::ARTIMTimer ( )`

## 27.14.3 Member Function Documentation

27.14.3.1 `double gdcm::network::ARTIMTimer::GetElapsedTime ( ) const`

27.14.3.2 `bool gdcm::network::ARTIMTimer::GetHasExpired ( ) const`

27.14.3.3 `double gdcm::network::ARTIMTimer::GetTimeout ( ) const`

27.14.3.4 `void gdcm::network::ARTIMTimer::SetTimeout ( double inTimeout )`

27.14.3.5 `void gdcm::network::ARTIMTimer::Start ( )`

27.14.3.6 `void gdcm::network::ARTIMTimer::Stop ( )`

The documentation for this class was generated from the following file:

- [gdcmARTIMTimer.h](#)

## 27.15 gdcm::ASN1 Class Reference

Class for [ASN1](#).

```
#include <gdcmASN1.h>
```

### Public Member Functions

- [ASN1](#) ()
- [~ASN1](#) ()

### Static Public Member Functions

- static bool [ParseDump](#) (const char \*array, size\_t length)
- static bool [ParseDumpFile](#) (const char \*filename)

### Protected Member Functions

- int [TestPBKDF2](#) ()

### 27.15.1 Detailed Description

Class for [ASN1](#).

### 27.15.2 Constructor & Destructor Documentation

27.15.2.1 `gdcm::ASN1::ASN1 ( )`

27.15.2.2 `gdcm::ASN1::~~ASN1 ( )`

### 27.15.3 Member Function Documentation

27.15.3.1 `static bool gdcm::ASN1::ParseDump ( const char * array, size_t length )` `[static]`

27.15.3.2 `static bool gdcm::ASN1::ParseDumpFile ( const char * filename )` `[static]`

27.15.3.3 `int gdcm::ASN1::TestPBKDF2 ( )` `[protected]`

The documentation for this class was generated from the following file:

- [gdcmASN1.h](#)

## 27.16 gdcm::network::AsynchronousOperationsWindowSub Class Reference

[AsynchronousOperationsWindowSub](#) PS 3.7 [Table D.3-7 ASYNCHRONOUS OPERATIONS WINDOW SUB-ITEM F](#)↵  
IELDS (A-ASSOCIATE-RQ)

```
#include <gdcmAsynchronousOperationsWindowSub.h>
```

### Public Member Functions

- [AsynchronousOperationsWindowSub](#) ( )
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) ( ) const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.16.1 Detailed Description

[AsynchronousOperationsWindowSub](#) PS 3.7 [Table D.3-7 ASYNCHRONOUS OPERATIONS WINDOW SUB-ITEM F](#)↵  
IELDS (A-ASSOCIATE-RQ)

### 27.16.2 Constructor & Destructor Documentation

27.16.2.1 `gdcm::network::AsynchronousOperationsWindowSub::AsynchronousOperationsWindowSub ( )`

### 27.16.3 Member Function Documentation

27.16.3.1 void gdcm::network::AsynchronousOperationsWindowSub::Print ( std::ostream & os ) const

27.16.3.2 std::istream& gdcm::network::AsynchronousOperationsWindowSub::Read ( std::istream & is )

27.16.3.3 size\_t gdcm::network::AsynchronousOperationsWindowSub::Size ( ) const

27.16.3.4 const std::ostream& gdcm::network::AsynchronousOperationsWindowSub::Write ( std::ostream & os ) const

The documentation for this class was generated from the following file:


- [gdcmAsynchronousOperationsWindowSub.h](#)

## 27.17 gdcm::Attribute< Group, Element, TVR, TVM > Class Template Reference

[Attribute](#) class This class use template metaprograming tricks to let the user know when the template instantiation does not match the public dictionary.

```
#include <gdcmAttribute.h>
```

Collaboration diagram for gdcm::Attribute< Group, Element, TVR, TVM >:

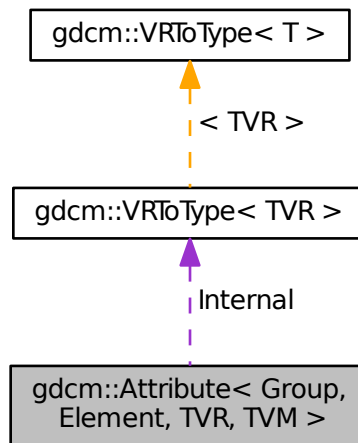

### Public Types

- enum { `VMType` = `VMToLength<TVM>::Length` }
- typedef `VRTToType< TVR >::Type` `ArrayType`

### Public Member Functions

- `GDCM_STATIC_ASSERT` (((`VR::VRType`) TVR &(`VR::VRType`)(`TagToType< Group, Element >::VRType`))))

- `GDCM_STATIC_ASSERT` (((`VM::VMType`) TVM &(`VM::VMType`)(TagToType< Group, `Element` >::`VMType`)))
- `GDCM_STATIC_ASSERT` (((((`VR::VRType`) TVR &`VR::VR_VM1`)&&((`VM::VMType`) TVM==`VM::VM1`))||!((`VR::VRType`) TVR &`VR::VR_VM1`)))
- `DataElement` `GetAsDataElement` () const
- unsigned int `GetNumberOfValues` () const
- `ArrayType` & `GetValue` (unsigned int idx=0)
- `ArrayType` const & `GetValue` (unsigned int idx=0) const
- const `ArrayType` \* `GetValues` () const
- bool `operator!=` (const `Attribute` &att) const
- bool `operator<` (const `Attribute` &att) const
- bool `operator==` (const `Attribute` &att) const
- `ArrayType` & `operator[]` (unsigned int idx)
- `ArrayType` const & `operator[]` (unsigned int idx) const
- void `Print` (std::ostream &os) const
- void `Set` (`DataSet` const &ds)
- void `SetFromDataElement` (`DataElement` const &de)
- void `SetFromDataSet` (`DataSet` const &ds)
- void `SetValue` (`ArrayType` v, unsigned int idx=0)
- void `SetValues` (const `ArrayType` \*array, unsigned int numel=`VMType`)

### Static Public Member Functions

- static `VM` `GetDictVM` ()
- static `VR` `GetDictVR` ()
- static `Tag` `GetTag` ()
- static `VM` `GetVM` ()
- static `VR` `GetVR` ()

### Public Attributes

- `ArrayType` `Internal` [`VMToLength`< `TVM` >::Length]

### Protected Member Functions

- void `SetByteValue` (const `ByteValue` \*bv)
- void `SetByteValueNoSwap` (const `ByteValue` \*bv)

## 27.17.1 Detailed Description

`template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType>class gdcmm::Attribute< Group, Element, TVR, TVM >`

`Attribute` class This class use template metaprograming tricks to let the user know when the template instantiation does not match the public dictionary.

Typical example that compile is: `Attribute<0x0008,0x9007> a = {"ORIGINAL","PRIMARY","T1","NONE"};`

Examples that will NOT compile are:

`Attribute<0x0018,0x1182, VR::IS, VM::VM1> fd1 = {};` // not enough parameters `Attribute<0x0018,0x1182, VR::IS, VM::VM2> fd2 = {0,1,2};` // too many initializers `Attribute<0x0018,0x1182, VR::IS, VM::VM3> fd3 = {0,1,2};` // VM3 is not valid `Attribute<0x0018,0x1182, VR::UL, VM::VM2> fd3 = {0,1};` // UL is not valid `VR`

Examples:

[CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream←\\_Image\\_Writer.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenFakeIdentifyFile.cxx](#), [GetSequenceUltrasound.cxx](#), [HelloWorld.cxx](#), [LargeVRDSExplicit.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadAndPrintAttributes.cxx](#), [rle2img.cxx](#), [SortImage.cxx](#), [StreamImageReaderTest.cxx](#), and [VolumeSorter.cxx](#).

## 27.17.2 Member Typedef Documentation

27.17.2.1 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> typedef VRToType<TVR>::Type gdcm::Attribute< Group, Element, TVR, TVM >::ArrayType`

## 27.17.3 Member Enumeration Documentation

27.17.3.1 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> anonymous enum`

Enumerator

***VMType***

## 27.17.4 Member Function Documentation

27.17.4.1 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> gdcm::Attribute< Group, Element, TVR, TVM >::GDCM_STATIC_ASSERT ( ((VR::VRType) TVR & (VR::VRType)(TagToType< Group, Element >::VRType)) )`

27.17.4.2 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> gdcm::Attribute< Group, Element, TVR, TVM >::GDCM_STATIC_ASSERT ( ((VM::VMType) TVM & (VM::VMType)(TagToType< Group, Element >::VMType)) )`

27.17.4.3 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> gdcm::Attribute< Group, Element, TVR, TVM >::GDCM_STATIC_ASSERT ( (((VR::VRType) TVR & VR::VR_VM1)&&((VM::VMType) TVM==VM::VM1))||!((VR::VRType) TVR & VR::VR_VM1)) )`

27.17.4.4 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> DataElement gdcm::Attribute< Group, Element, TVR, TVM >::GetAsDataElement ( ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcm::DataElement::GetVR()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetVR()`, `gdcm::DataElement::SetByteValue()`, `gdcm::DataElement::SetVR()`, `gdcm::VR::SQ`, `gdcm::VR::UI`, and `gdcm::VR::V←RASCII`.

27.17.4.5 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> static VM gdcm::Attribute< Group, Element, TVR, TVM >::GetDictVM ( ) [inline], [static]`

27.17.4.6 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> static VR gdcmm::Attribute< Group, Element, TVR, TVM >::GetDictVR ( ) [inline], [static]`

27.17.4.7 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> unsigned int gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues ( ) const [inline]`

Referenced by `gdcmm::Attribute< Group, Element, TVR, TVM >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetValue()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GetValue()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::operator!=()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::operator!=()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::operator<()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::operator<()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::operator==()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::operator==()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::Print()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::Print()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetByteValue()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValue()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValueNoSwap()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetValue()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetValue()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetValues()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetValues()`.

27.17.4.8 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> static Tag gdcmm::Attribute< Group, Element, TVR, TVM >::GetTag ( ) [inline], [static]`

Referenced by `gdcmm::Attribute< Group, Element, TVR, TVM >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::Print()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::Print()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::Print()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::Set()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::Set()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::Set()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataSet()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataSet()`.

27.17.4.9 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> ArrayType& gdcmm::Attribute< Group, Element, TVR, TVM >::GetValue ( unsigned int idx = 0 ) [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`.

Referenced by `gdcmm::Attribute< Group, Element, TVR, TVM >::operator[]()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::operator[]()`.

27.17.4.10 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> ArrayType const& gdcmm::Attribute< Group, Element, TVR, TVM >::GetValue ( unsigned int idx = 0 ) const [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`.



```
27.17.4.11 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> const ArrayType* gdcmm::Attribute< Group, Element, TVR, TVM
>::GetValues ( ) const [inline]
```

References gdcmm::Attribute< Group, Element, TVR, TVM >::Internal.

Referenced by gdcmm::Attribute< Group, Element, TVR, TVM >::operator!==( ), gdcmm::Attribute< Group, Element, T←VR, VM::VM1 >::operator!==( ), gdcmm::Attribute< Group, Element, TVR, TVM >::operator<(), gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::operator<(), gdcmm::Attribute< Group, Element, TVR, TVM >::operator==( ), and gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::operator==( ).

```
27.17.4.12 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> static VM gdcmm::Attribute< Group, Element, TVR, TVM >::GetVM ( )
[inline], [static]
```

Referenced by gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >::GetDictVM(), and gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >::Print().

```
27.17.4.13 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> static VR gdcmm::Attribute< Group, Element, TVR, TVM >::GetVR ( )
[inline], [static]
```

Referenced by gdcmm::Attribute< Group, Element, TVR, TVM >::GetAsDataElement(), gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement(), gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >::GetAsData←Element(), gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >::Print(), gdcmm::Attribute< Group, Element, TV←R, TVM >::SetFromDataElement(), gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement(), and gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >::SetFromDataElement().

```
27.17.4.14 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> bool gdcmm::Attribute< Group, Element, TVR, TVM >::operator!= ( const
Attribute< Group, Element, TVR, TVM > & att ) const [inline]
```

References gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues(), and gdcmm::Attribute< Group, Element, TVR, TVM >::GetValues().

```
27.17.4.15 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> bool gdcmm::Attribute< Group, Element, TVR, TVM >::operator< ( const
Attribute< Group, Element, TVR, TVM > & att ) const [inline]
```

References gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues(), and gdcmm::Attribute< Group, Element, TVR, TVM >::GetValues().

```
27.17.4.16 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> bool gdcmm::Attribute< Group, Element, TVR, TVM >::operator==( const
Attribute< Group, Element, TVR, TVM > & att ) const [inline]
```

References gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues(), and gdcmm::Attribute< Group, Element, TVR, TVM >::GetValues().

27.17.4.17 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> ArrayType& gdcmm::Attribute< Group, Element, TVR, TVM >::operator[] ( unsigned int idx ) [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetValue()`.

27.17.4.18 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> ArrayType const& gdcmm::Attribute< Group, Element, TVR, TVM >::operator[] ( unsigned int idx ) const [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetValue()`.

27.17.4.19 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> void gdcmm::Attribute< Group, Element, TVR, TVM >::Print ( std::ostream & os ) const [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcmm::Attribute< Group, Element, TVR, TVM >::GetTag()`.

27.17.4.20 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> void gdcmm::Attribute< Group, Element, TVR, TVM >::Set ( DataSet const & ds ) [inline]`

References `gdcmm::DataSet::GetDataElement()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetTag()`, and `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`.

27.17.4.21 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> void gdcmm::Attribute< Group, Element, TVR, TVM >::SetByteValue ( const ByteValue * bv ) [inline], [protected]`

References `gdcmm::ByteValue::GetLength()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcmm::ByteValue::GetPointer()`.

Referenced by `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`.

27.17.4.22 `template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM = TagToType<Group, Element>::VMType> void gdcmm::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap ( const ByteValue * bv ) [inline], [protected]`

References `gdcmm::ByteValue::GetLength()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcmm::ByteValue::GetPointer()`.

Referenced by `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`.

```
27.17.4.23 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM
= TagToType<Group, Element>::VMType> void gdcm::Attribute< Group, Element, TVR, TVM
>::SetFromDataElement ( DataElement const & de ) [inline]
```

References gdcm::DataElement::GetByteValue(), gdcm::Tag::GetGroup(), gdcm::DataElement::GetTag(), gdcm::Attribute< Group, Element, TVR, TVM >::GetTag(), gdcm::DataElement::GetVR(), gdcm::Attribute< Group, Element, TVR, TVM >::GetVR(), gdcm::VR::INVALID, gdcm::DataElement::IsEmpty(), gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValue(), gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap(), and gdcm::VR::UN.

Referenced by gdcm::Attribute< Group, Element, TVR, TVM >::Set(), gdcm::Attribute< Group, Element, TVR, VM::VM1 >::Set(), gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::Set(), gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataSet(), gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet(), and gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::SetFromDataSet().

```
27.17.4.24 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> void gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataSet (
DataSet const & ds ) [inline]
```

References gdcm::DataSet::FindDataElement(), gdcm::DataSet::GetDataElement(), gdcm::Attribute< Group, Element, TVR, TVM >::GetTag(), gdcm::DataElement::IsEmpty(), and gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement().

```
27.17.4.25 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> void gdcm::Attribute< Group, Element, TVR, TVM >::SetValue (
ArrayType v, unsigned int idx = 0 ) [inline]
```

References gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues().

```
27.17.4.26 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> void gdcm::Attribute< Group, Element, TVR, TVM >::SetValues ( const
ArrayType * array, unsigned int numel = VMType ) [inline]
```

Examples:

[LargeVRDSExplicit.cxx](#).

References gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues().

Referenced by gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::SetByteValue(), and gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::SetNumberOfValues().

## 27.17.5 Member Data Documentation

```
27.17.5.1 template<uint16_t Group, uint16_t Element, int TVR = TagToType<Group, Element>::VRType, int TVM =
TagToType<Group, Element>::VMType> ArrayType gdcm::Attribute< Group, Element, TVR, TVM
>::Internal[VMToLength< TVM >::Length]
```

Referenced by gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetValue(), gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetValues(), gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::GetValues(), gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::Print(), gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::SetValues(), and gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >::~~Attribute().


The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)

## 27.18 gdcm::Attribute< Group, Element, TVR, VM::VM1 > Class Template Reference

```
#include <gdcmAttribute.h>
```

Collaboration diagram for gdcm::Attribute< Group, Element, TVR, VM::VM1 >:

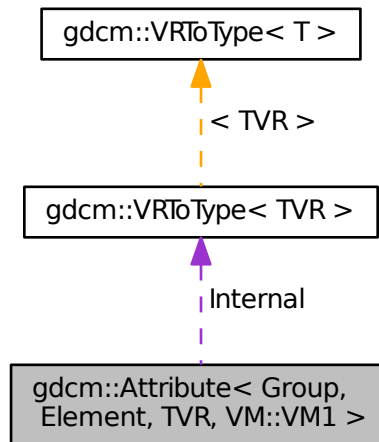

### Public Types

- enum { `VMType` = `VMToLength<VM::VM1>::Length` }
- typedef `VRTToType< TVR >::Type` `ArrayType`

### Public Member Functions

- `GDCM_STATIC_ASSERT (VMToLength< VM::VM1 >::Length==1)`
- `GDCM_STATIC_ASSERT (((VR::VRTType) TVR &(VR::VRTType)(TagToType< Group, Element >::VRTType)))`
- `GDCM_STATIC_ASSERT (((VM::VMType) VM::VM1 &(VM::VMType)(TagToType< Group, Element >::VMType)))`
- `GDCM_STATIC_ASSERT (((((VR::VRTType) TVR &VR::VR_VM1)&((VM::VMType) VM::VM1==VM::VM1))||((VR::VRTType) TVR &VR::VR_VM1)))`
- `DataElement GetAsDataElement () const`
- `unsigned int GetNumberOfValues () const`
- `ArrayType & GetValue ()`
- `ArrayType const & GetValue () const`
- `const ArrayType * GetValues () const`

- bool [operator!=](#) (const [Attribute](#) &att) const
- bool [operator<](#) (const [Attribute](#) &att) const
- bool [operator==](#) (const [Attribute](#) &att) const
- void [Print](#) (std::ostream &os) const
- void [Set](#) ([DataSet](#) const &ds)
- void [SetFromDataElement](#) ([DataElement](#) const &de)
- void [SetFromDataSet](#) ([DataSet](#) const &ds)
- void [SetValue](#) ([ArrayType](#) v)

### Static Public Member Functions

- static [VM GetDictVM](#) ()
- static [VR GetDictVR](#) ()
- static [Tag GetTag](#) ()
- static [VM GetVM](#) ()
- static [VR GetVR](#) ()

### Public Attributes

- [ArrayType Internal](#)

### Protected Member Functions

- void [SetByteValue](#) (const [ByteValue](#) \*bv)
- void [SetByteValueNoSwap](#) (const [ByteValue](#) \*bv)

## 27.18.1 Member Typedef Documentation

- 27.18.1.1 `template<uint16_t Group, uint16_t Element, int TVR> typedef VRToType<TVR>::Type gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::ArrayType`

## 27.18.2 Member Enumeration Documentation

- 27.18.2.1 `template<uint16_t Group, uint16_t Element, int TVR> anonymous enum`

Enumerator

***VMType***

## 27.18.3 Member Function Documentation

- 27.18.3.1 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GDCM_STATIC_ASSERT ( VMToLength< VM::VM1 >::Length == 1 )`
- 27.18.3.2 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GDCM_STATIC_ASSERT ( ((VR::VRType) TVR &(VR::VRType)(TagToType< Group, Element >::VRType)) )`

27.18.3.3 `template<uint16_t Group, uint16_t Element, int TVR> gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GDCM_STATIC_ASSERT ( ((VM::VMType) VM::VM1 &(VM::VMType)(TagToType< Group, Element >::VMType)) )`

27.18.3.4 `template<uint16_t Group, uint16_t Element, int TVR> gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GDCM_STATIC_ASSERT ( (((VR::VRType) TVR &VR::VR_VM1)&&((VM::VMType) VM::VM1==VM::VM1))||!((VR::VRType) TVR &VR::VR_VM1)) )`

27.18.3.5 `template<uint16_t Group, uint16_t Element, int TVR> DataElement gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement ( ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcm::DataElement::GetVR()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetVR()`, `gdcm::DataElement::SetByteValue()`, `gdcm::DataElement::SetVR()`, `gdcm::VR::SQ`, `gdcm::VR::UI`, and `gdcm::VR::VRASCII`.

27.18.3.6 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetDictVM ( ) [inline],[static]`

27.18.3.7 `template<uint16_t Group, uint16_t Element, int TVR> static VR gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetDictVR ( ) [inline],[static]`

27.18.3.8 `template<uint16_t Group, uint16_t Element, int TVR> unsigned int gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetNumberOfValues ( ) const [inline]`

27.18.3.9 `template<uint16_t Group, uint16_t Element, int TVR> static Tag gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetTag ( ) [inline],[static]`

27.18.3.10 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType& gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetValue ( ) [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::Internal`.

27.18.3.11 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType const& gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetValue ( ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::Internal`.

27.18.3.12 `template<uint16_t Group, uint16_t Element, int TVR> const ArrayType* gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetValues ( ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::Internal`.

27.18.3.13 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetVM ( ) [inline],[static]`

References `gdcm::VM::VM1`.

27.18.3.14 `template<uint16_t Group, uint16_t Element, int TVR> static VR gdcm::Attribute< Group, Element, TVR, VM::VM1 >::GetVR( ) [inline], [static]`

27.18.3.15 `template<uint16_t Group, uint16_t Element, int TVR> bool gdcm::Attribute< Group, Element, TVR, VM::VM1 >::operator!=( const Attribute< Group, Element, TVR, VM::VM1 > & att ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcm::Attribute< Group, Element, TVR, TVM >::GetValues()`.

27.18.3.16 `template<uint16_t Group, uint16_t Element, int TVR> bool gdcm::Attribute< Group, Element, TVR, VM::VM1 >::operator< ( const Attribute< Group, Element, TVR, VM::VM1 > & att ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcm::Attribute< Group, Element, TVR, TVM >::GetValues()`.

27.18.3.17 `template<uint16_t Group, uint16_t Element, int TVR> bool gdcm::Attribute< Group, Element, TVR, VM::VM1 >::operator==( const Attribute< Group, Element, TVR, VM::VM1 > & att ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcm::Attribute< Group, Element, TVR, TVM >::GetValues()`.

27.18.3.18 `template<uint16_t Group, uint16_t Element, int TVR> void gdcm::Attribute< Group, Element, TVR, VM::VM1 >::Print( std::ostream & os ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetTag()`, and `gdcm::Attribute< Group, Element, TVR, TVM >::Internal`.

27.18.3.19 `template<uint16_t Group, uint16_t Element, int TVR> void gdcm::Attribute< Group, Element, TVR, VM::VM1 >::Set( DataSet const & ds ) [inline]`

References `gdcm::DataSet::GetDataElement()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetTag()`, and `gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`.

27.18.3.20 `template<uint16_t Group, uint16_t Element, int TVR> void gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValue( const ByteValue * bv ) [inline], [protected]`

References `gdcm::ByteValue::GetLength()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcm::ByteValue::GetPointer()`.

27.18.3.21 `template<uint16_t Group, uint16_t Element, int TVR> void gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValueNoSwap( const ByteValue * bv ) [inline], [protected]`

References `gdcm::ByteValue::GetLength()`, `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, and `gdcm::ByteValue::GetPointer()`.

27.18.3.22 `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement ( DataElement const & de ) [inline]`

References `gdcM::DataElement::GetByteValue()`, `gdcM::Tag::GetGroup()`, `gdcM::DataElement::GetTag()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcM::DataElement::GetVR()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetVR()`, `gdcM::VR::INVALID`, `gdcM::DataElement::IsEmpty()`, `gdcM::Attribute< Group, Element, TVR, TVM >::SetByteValue()`, `gdcM::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap()`, and `gdcM::VR::UN`.

27.18.3.23 `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet ( DataSet const & ds ) [inline]`

References `gdcM::DataSet::FindDataElement()`, `gdcM::DataSet::GetDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcM::DataElement::IsEmpty()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`.

27.18.3.24 `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetValue ( ArrayType v ) [inline]`

## 27.18.4 Member Data Documentation

27.18.4.1 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType gdcM::Attribute< Group, Element, TVR, VM::VM1 >::Internal`


The documentation for this class was generated from the following file:

- [gdcMAttribute.h](#)

## 27.19 gdcM::Attribute< Group, Element, TVR, VM::VM1\_3 > Class Template Reference


```
#include <gdcMAttribute.h>
```

Inheritance diagram for `gdcM::Attribute< Group, Element, TVR, VM::VM1_3 >`:

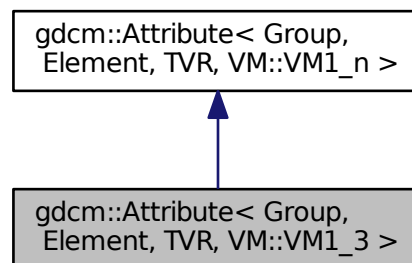



Collaboration diagram for `gdcm::Attribute< Group, Element, TVR, VM::VM1_3 >`:

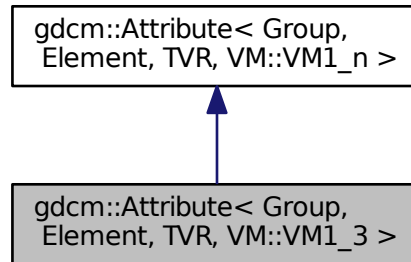

## Public Member Functions

- [VM GetVM](#) () const

## Additional Inherited Members

### 27.19.1 Member Function Documentation

27.19.1.1 `template<uint16_t Group, uint16_t Element, int TVR> VM gdcm::Attribute< Group, Element, TVR, VM::VM1_3 >::GetVM ( ) const` [inline]

References `gdcm::VM::VM1_3`.


The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)


## 27.20 `gdcm::Attribute< Group, Element, TVR, VM::VM1_8 >` Class Template Reference

```
#include <gdcmAttribute.h>
```

Inheritance diagram for `gdc::Attribute< Group, Element, TVR, VM::VM1_8 >`:

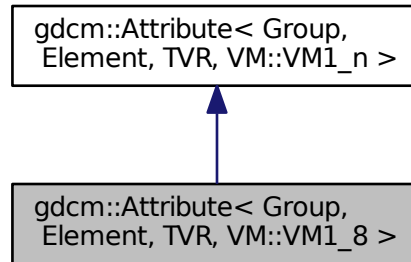

Collaboration diagram for `gdc::Attribute< Group, Element, TVR, VM::VM1_8 >`:

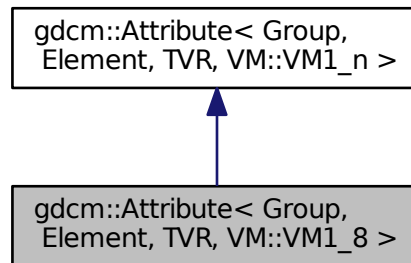

## Public Member Functions

- [VM GetVM](#) () const

## Additional Inherited Members

### 27.20.1 Member Function Documentation

27.20.1.1 `template<uint16_t Group, uint16_t Element, int TVR> VM gdc::Attribute< Group, Element, TVR, VM::VM1_8 >::GetVM ( ) const [inline]`

References `gdc::VM::VM1_8`.


The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)

## 27.21 gdcm::Attribute< Group, Element, TVR, VM::VM1\_n > Class Template Reference

```
#include <gdcmAttribute.h>
```

Inheritance diagram for gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >:

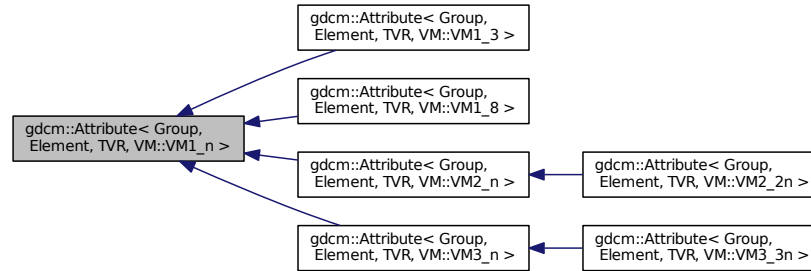

### Public Types

- typedef `VRToType< TVR >::Type` `ArrayType`

### Public Member Functions

- `Attribute ()`
- `~Attribute ()`
- `GDCM_STATIC_ASSERT (((VR::VRType) TVR & (VR::VRType)(TagToType< Group, Element >::VRType)))`
- `GDCM_STATIC_ASSERT ((VM::VM1_n & (VM::VMType)(TagToType< Group, Element >::VMType)))`
- `GDCM_STATIC_ASSERT (((((VR::VRType) TVR & VR::VR_VM1) && ((VM::VMType) TagToType< Group, Element >::VMType==VM::VM1)) || !((VR::VRType) TVR & VR::VR_VM1)))`
- `DataElement GetAsDataElement () const`
- `unsigned int GetNumberOfValues () const`
- `ArrayType & GetValue (unsigned int idx=0)`
- `ArrayType const & GetValue (unsigned int idx=0) const`
- `const ArrayType * GetValues () const`
- `ArrayType & operator[] (unsigned int idx)`
- `ArrayType const & operator[] (unsigned int idx) const`
- `void Print (std::ostream &os) const`
- `void Set (DataSet const &ds)`
- `void SetFromDataElement (DataElement const &de)`
- `void SetFromDataSet (DataSet const &ds)`
- `void SetNumberOfValues (unsigned int numel)`
- `void SetValue (unsigned int idx, ArrayType v)`
- `void SetValue (ArrayType v)`
- `void SetValues (const ArrayType *array, unsigned int numel, bool own=false)`

## Static Public Member Functions

- static [VM GetDictVM](#) ()
- static [VR GetDictVR](#) ()
- static [Tag GetTag](#) ()
- static [VM GetVM](#) ()
- static [VR GetVR](#) ()

## Protected Member Functions

- void [SetByteValue](#) (const [ByteValue](#) \*bv)

### 27.21.1 Member Typedef Documentation

27.21.1.1 `template<uint16_t Group, uint16_t Element, int TVR> typedef VRToType<TVR>::Type gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::ArrayType`

### 27.21.2 Constructor & Destructor Documentation

27.21.2.1 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::Attribute ( ) [inline],[explicit]`

27.21.2.2 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::~~Attribute ( ) [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::Internal`.

### 27.21.3 Member Function Documentation

27.21.3.1 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GDCM_STATIC_ASSERT ( ((VR::VRType) TVR &(VR::VRType)(TagToType< Group, Element >::VRType)) )`

27.21.3.2 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GDCM_STATIC_ASSERT ( (VM::VM1_n &(VM::VMType)(TagToType< Group, Element >::VMType)) )`

27.21.3.3 `template<uint16_t Group, uint16_t Element, int TVR> gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GDCM_STATIC_ASSERT ( (((VR::VRType) TVR &VR::VR_VM1)&&((VM::VMType) TagToType< Group, Element >::VMType==VM::VM1))||!((VR::VRType) TVR &VR::VR_VM1)) )`

27.21.3.4 `template<uint16_t Group, uint16_t Element, int TVR> DataElement gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement ( ) const [inline]`

References `gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcmm::DataElement::GetVR()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::GetVR()`, `gdcmm::DataElement::SetByteValue()`, `gdcmm::DataElement::SetVR()`, `gdcmm::VR::SQ`, `gdcmm::VR::UI`, and `gdcmm::VR::V←RASCII`.

27.21.3.5 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetDictVM( ) [inline], [static]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetVM()`.

27.21.3.6 `template<uint16_t Group, uint16_t Element, int TVR> static VR gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetDictVR( ) [inline], [static]`

27.21.3.7 `template<uint16_t Group, uint16_t Element, int TVR> unsigned int gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetNumberOfValues( ) const [inline]`

27.21.3.8 `template<uint16_t Group, uint16_t Element, int TVR> static Tag gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetTag( ) [inline], [static]`

27.21.3.9 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType& gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetValue( unsigned int idx = 0 ) [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`.

27.21.3.10 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType const& gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetValue( unsigned int idx = 0 ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`.

27.21.3.11 `template<uint16_t Group, uint16_t Element, int TVR> const ArrayType* gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetValues( ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::Internal`.

27.21.3.12 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetVM( ) [inline], [static]`

References `gdcm::VM::VM1_n`.

27.21.3.13 `template<uint16_t Group, uint16_t Element, int TVR> static VR gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetVR( ) [inline], [static]`

27.21.3.14 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType& gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::operator[]( unsigned int idx ) [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetValue()`.

27.21.3.15 `template<uint16_t Group, uint16_t Element, int TVR> ArrayType const& gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::operator[]( unsigned int idx ) const [inline]`

References `gdcm::Attribute< Group, Element, TVR, TVM >::GetValue()`.

**27.21.3.16** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::Print ( std::ostream & os ) const [inline]`

References `gdcM::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetVM()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::GetVR()`.

**27.21.3.17** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::Set ( DataSet const & ds ) [inline]`

References `gdcM::DataSet::GetDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`.

**27.21.3.18** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetByteValue ( const ByteValue * bv ) [inline], [protected]`

References `gdcM::ByteValue::GetLength()`, `gdcM::ByteValue::GetPointer()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::SetValues()`.

**27.21.3.19** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement ( DataElement const & de ) [inline]`

References `gdcM::DataElement::GetByteValue()`, `gdcM::Tag::GetGroup()`, `gdcM::DataElement::GetTag()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcM::DataElement::GetVR()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetVR()`, `gdcM::DataElement::IsEmpty()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::SetByteValue()`.

**27.21.3.20** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataSet ( DataSet const & ds ) [inline]`

References `gdcM::DataSet::FindDataElement()`, `gdcM::DataSet::GetDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetTag()`, `gdcM::DataElement::IsEmpty()`, and `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`.

**27.21.3.21** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetNumberOfValues ( unsigned int numel ) [inline]`

References `gdcM::Attribute< Group, Element, TVR, TVM >::SetValues()`.

**27.21.3.22** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetValue ( unsigned int idx, ArrayType v ) [inline]`

References `gdcM::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues()`.

**27.21.3.23** `template<uint16_t Group, uint16_t Element, int TVR> void gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetValue ( ArrayType v ) [inline]`

References `SetValue()`.

Referenced by SetValue().

```
27.21.3.24 template<uint16_t Group, uint16_t Element, int TVR> void gdcmm::Attribute< Group, Element, TVR, VM::VM1_n
>::SetValues ( const ArrayType * array, unsigned int numel, bool own = false ) [inline]
```

References gdcmm::Attribute< Group, Element, TVR, TVM >::GetNumberOfValues(), and gdcmm::Attribute< Group, Element, TVR, TVM >::Internal.


The documentation for this class was generated from the following file:

- [gdcmmAttribute.h](#)


## 27.22 gdcmm::Attribute< Group, Element, TVR, VM::VM2\_2n > Class Template Reference

```
#include <gdcmmAttribute.h>
```

Inheritance diagram for gdcmm::Attribute< Group, Element, TVR, VM::VM2\_2n >:

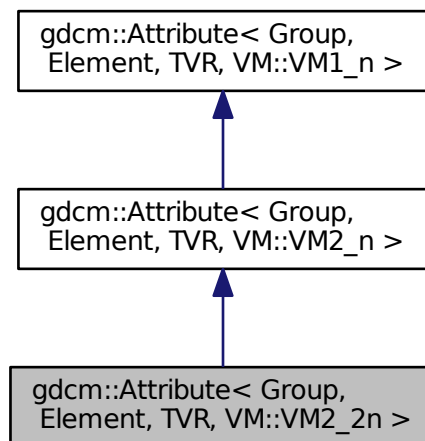

Collaboration diagram for `gdcM::Attribute< Group, Element, TVR, VM::VM2_2n >`:

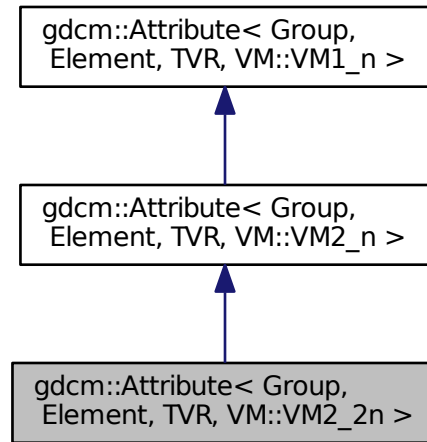

## Static Public Member Functions

- static [VM GetVM](#) ()

## Additional Inherited Members

### 27.22.1 Member Function Documentation

27.22.1.1 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcM::Attribute< Group, Element, TVR, VM::VM2_2n >::GetVM ( ) [inline], [static]`

References `gdcM::VM::VM2_2n`.

The documentation for this class was generated from the following file:


- [gdcMAttribute.h](#)

## 27.23 `gdcM::Attribute< Group, Element, TVR, VM::VM2_n >` Class Template Reference


```
#include <gdcMAttribute.h>
```



Inheritance diagram for gdcm::Attribute< Group, Element, TVR, VM::VM2\_n >:

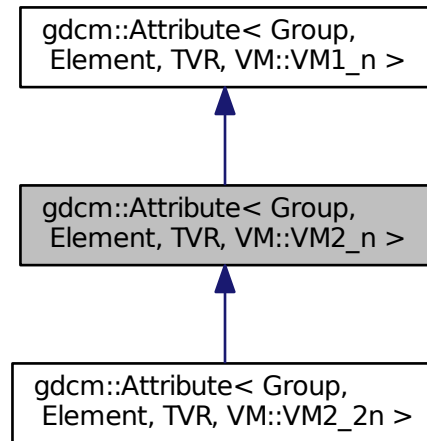

Collaboration diagram for gdcm::Attribute< Group, Element, TVR, VM::VM2\_n >:

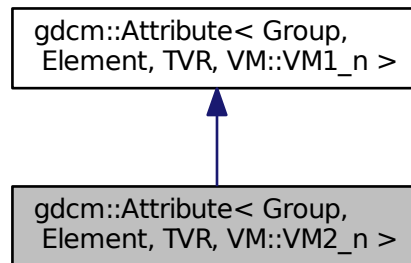

### Public Member Functions

- [VM GetVM](#) () const

### Additional Inherited Members

#### 27.23.1 Member Function Documentation

27.23.1.1 `template<uint16_t Group, uint16_t Element, int TVR> VM gdcM::Attribute< Group, Element, TVR, VM::VM2_n >::GetVM( ) const [inline]`

References `gdcM::VM::VM2_n`.


The documentation for this class was generated from the following file:

- [gdcMAttribute.h](#)


## 27.24 `gdcM::Attribute< Group, Element, TVR, VM::VM3_3n >` Class Template Reference

```
#include <gdcMAttribute.h>
```

Inheritance diagram for `gdcM::Attribute< Group, Element, TVR, VM::VM3_3n >`:

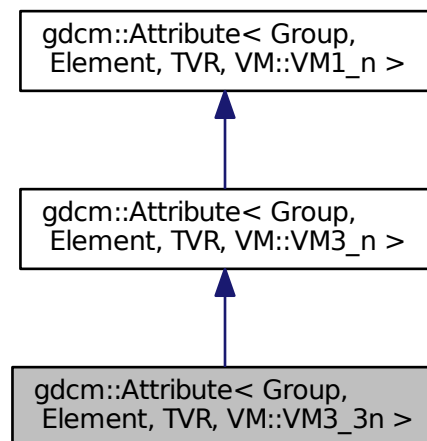

Collaboration diagram for gdcmm::Attribute< Group, Element, TVR, VM::VM3\_3n >:

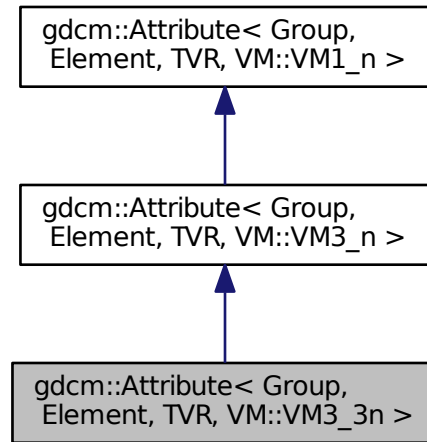

### Static Public Member Functions

- static [VM GetVM](#) ()

### Additional Inherited Members

#### 27.24.1 Member Function Documentation

27.24.1.1 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcmm::Attribute< Group, Element, TVR, VM::VM3_3n >::GetVM ( ) [inline], [static]`

References `gdcmm::VM::VM3_3n`.


The documentation for this class was generated from the following file:

- [gdcmmAttribute.h](#)


## 27.25 gdcmm::Attribute< Group, Element, TVR, VM::VM3\_n > Class Template Reference

```
#include <gdcmmAttribute.h>
```

Inheritance diagram for `gdcM::Attribute< Group, Element, TVR, VM::VM3_n >`:

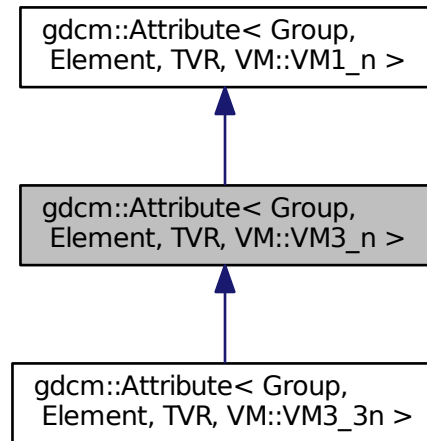

Collaboration diagram for `gdcM::Attribute< Group, Element, TVR, VM::VM3_n >`:

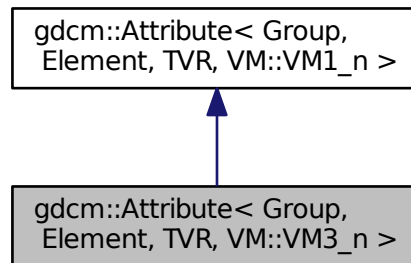

### Static Public Member Functions

- static [VM GetVM](#) ()

### Additional Inherited Members

#### 27.25.1 Member Function Documentation

27.25.1.1 `template<uint16_t Group, uint16_t Element, int TVR> static VM gdcm::Attribute< Group, Element, TVR, VM::VM3_n >::GetVM( ) [inline],[static]`

References `gdcm::VM::VM3_n`.

The documentation for this class was generated from the following file:


- [gdcmAttribute.h](#)

## 27.26 gdcm::AudioCodec Class Reference


[AudioCodec](#).

```
#include <gdcmAudioCodec.h>
```

Inheritance diagram for `gdcm::AudioCodec`:

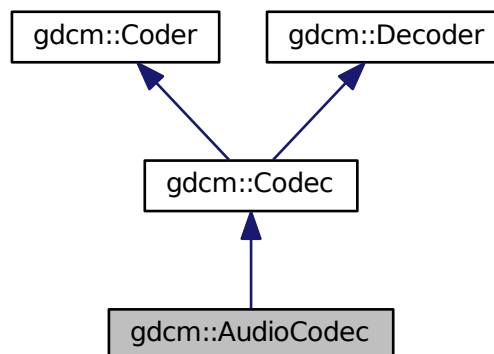

Collaboration diagram for `gdcm::AudioCodec`:

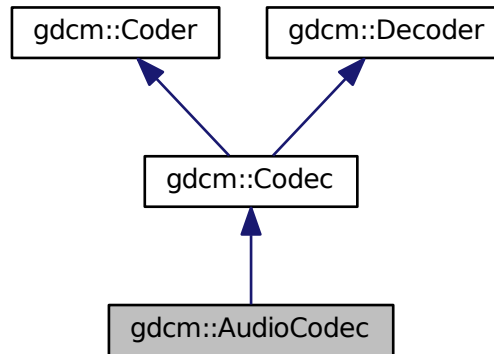

## Public Member Functions

- [AudioCodec](#) ()
- [~AudioCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &) const  
*Return whether this coder support this transfer syntax (can code it)*
- bool [CanDecode](#) ([TransferSyntax](#) const &) const  
*Return whether this decoder support this transfer syntax (can decode it)*
- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)  
*Decode.*

## Additional Inherited Members

### 27.26.1 Detailed Description

[AudioCodec](#).

### 27.26.2 Constructor & Destructor Documentation

27.26.2.1 `gdcm::AudioCodec::AudioCodec ( )`

27.26.2.2 `gdcm::AudioCodec::~~AudioCodec ( )`

### 27.26.3 Member Function Documentation

27.26.3.1 `bool gdcm::AudioCodec::CanCode ( TransferSyntax const & ) const` `[inline]`, `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Implements [gdcm::Coder](#).

27.26.3.2 `bool gdcm::AudioCodec::CanDecode ( TransferSyntax const & ) const` `[inline], [virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Implements [gdcm::Decoder](#).

27.26.3.3 `bool gdcm::AudioCodec::Decode ( DataElement const & , DataElement & )` `[virtual]`

Decode.

Reimplemented from [gdcm::Decoder](#).

The documentation for this class was generated from the following file:

- [gdcmAudioCodec.h](#)

## 27.27 gdcm::Base64 Class Reference

Class for [Base64](#).

```
#include <gdcmBase64.h>
```

### Static Public Member Functions

- static `size_t` [Decode](#) (`char *dst`, `size_t dlen`, `const char *src`, `size_t slen`)  
*Decode a base64-formatted buffer.*
- static `size_t` [Encode](#) (`char *dst`, `size_t dlen`, `const char *src`, `size_t slen`)  
*Encode a buffer into base64 format.*
- static `size_t` [GetDecodeLength](#) (`const char *src`, `size_t len`)
- static `size_t` [GetEncodeLength](#) (`const char *src`, `size_t srclen`)

### 27.27.1 Detailed Description

Class for [Base64](#).

### 27.27.2 Member Function Documentation

27.27.2.1 `static size_t gdcm::Base64::Decode ( char * dst, size_t dlen, const char * src, size_t slen )` `[static]`

Decode a base64-formatted buffer.

Parameters

<code>dst</code>	destination buffer
------------------	--------------------

<i>dlen</i>	size of the buffer
<i>src</i>	source buffer
<i>slen</i>	amount of data to be decoded

**Returns**

0 if not successful, size of decoded otherwise

**Examples:**

[DumpExamCard.cxx](#).

**27.27.2.2** `static size_t gdcm::Base64::Encode ( char * dst, size_t dlen, const char * src, size_t slen ) [static]`

Encode a buffer into base64 format.

**Parameters**

<i>dst</i>	destination buffer
<i>dlen</i>	size of the buffer
<i>src</i>	source buffer
<i>slen</i>	amount of data to be encoded

**Returns**

0 if not successful, size of encoded otherwise

**27.27.2.3** `static size_t gdcm::Base64::GetDecodeLength ( const char * src, size_t len ) [static]`

Call this function to obtain the required buffer size

**Examples:**

[DumpExamCard.cxx](#).

**27.27.2.4** `static size_t gdcm::Base64::GetEncodeLength ( const char * src, size_t srclen ) [static]`

Call this function to obtain the required buffer size

The documentation for this class was generated from the following file:

- [gdcmBase64.h](#)


## 27.28 gdcm::network::BaseCompositeMessage Class Reference

[BaseCompositeMessage](#) The Composite events described in section 3.7-2009 of the DICOM standard all use their own messages. These messages are constructed using Presentation Data Values, from section 3.8-2009 of the standard, and then fill in appropriate values in their datasets.

```
#include <gdcmBaseCompositeMessage.h>
```



Inheritance diagram for gdcmm::network::BaseCompositeMessage:

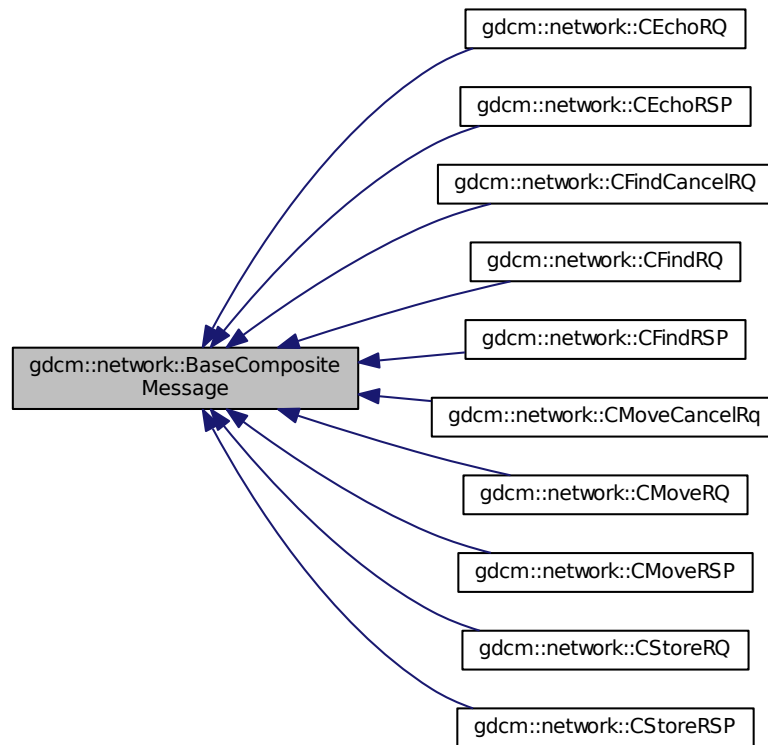

### Public Member Functions

- virtual `~BaseCompositeMessage()`
- virtual `std::vector< PresentationDataValue > ConstructPDV (const ULConnection &inConnection, const Base←RootQuery *inRootQuery)=0`

### 27.28.1 Detailed Description

**BaseCompositeMessage** The Composite events described in section 3.7-2009 of the DICOM standard all use their own messages. These messages are constructed using Presentation Data Values, from section 3.8-2009 of the standard, and then fill in appropriate values in their datasets.

So, for the five composites:

- C-ECHO
- C-FIND
- C-MOVE
- C-GET

- C-STORE there are a series of messages. However, all of these messages are obtained as part of a PDataPDU, and all have to be placed there. Therefore, since they all have shared functionality and construction tropes, that will be put into a base class. Further, the base class will be then returned by the factory class, `gdcmCompositePDUFactory`.

This is an abstract class. It cannot be instantiated on its own.

## 27.28.2 Constructor & Destructor Documentation

27.28.2.1 `virtual gdcm::network::BaseCompositeMessage::~BaseCompositeMessage ( ) [inline], [virtual]`

## 27.28.3 Member Function Documentation

27.28.3.1 `virtual std::vector<PresentationDataValue> gdcm::network::BaseCompositeMessage::ConstructPDV ( const ULConnection & inConnection, const BaseRootQuery * inRootQuery ) [pure virtual]`

Implemented in [gdcm::network::CMoveRQ](#), [gdcm::network::CFindRQ](#), and [gdcm::network::CEchoRQ](#).

The documentation for this class was generated from the following file:


- [gdcmBaseCompositeMessage.h](#)

## 27.29 gdcm::network::BasePDU Class Reference

[BasePDU](#) base class for PDUs.

```
#include <gdcmBasePDU.h>
```

Inheritance diagram for gdcmm::network::BasePDU:

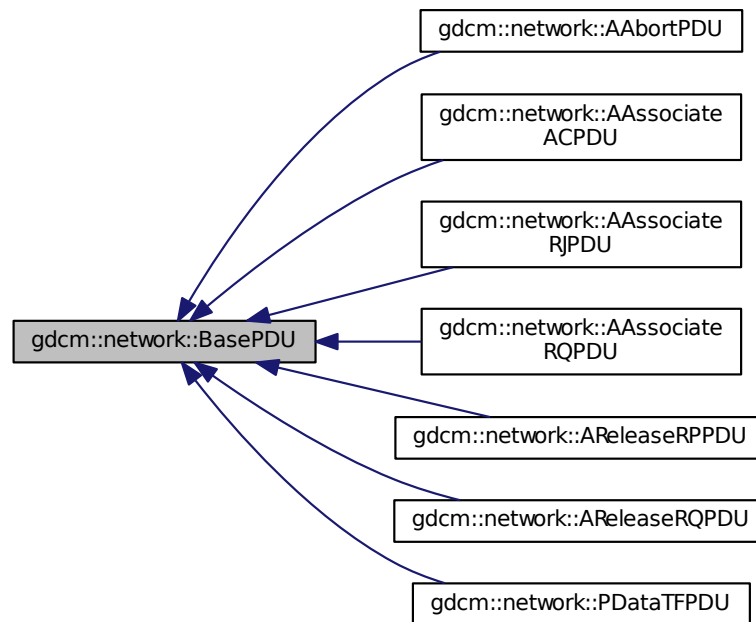

## Public Member Functions

- virtual [~BasePDU](#) ()
- virtual bool [IsLastFragment](#) () const =0
- virtual void [Print](#) (std::ostream &os) const =0
- virtual std::istream & [Read](#) (std::istream &is)=0
- virtual size\_t [Size](#) () const =0
- virtual const std::ostream & [Write](#) (std::ostream &os) const =0

### 27.29.1 Detailed Description

[BasePDU](#) base class for PDUs.

all PDUs start with the first ten bytes as specified: 01 PDU type 02 reserved 3-6 PDU Length (unsigned) 7-10 variable on some, 7-10 are split (7-8 as protocol version in Associate-RQ, for instance, while associate-rj splits those four bytes differently).

Also common to all the PDUs is their ability to read and write to a stream.

So, let's just get them all bunched together into one (abstract) class, shall we?

Why? 1) so that the [ULEvent](#) can have the PDU stored in it, since the event takes PDUs and not other class structures (other class structures get converted into PDUs) 2) to make reading PDUs in the event loop cleaner

## 27.29.2 Constructor & Destructor Documentation

27.29.2.1 `virtual gdcmm::network::BasePDU::~BasePDU ( ) [inline], [virtual]`

## 27.29.3 Member Function Documentation

27.29.3.1 `virtual bool gdcmm::network::BasePDU::IsLastFragment ( ) const [pure virtual]`

Implemented in [gdcmm::network::AAssociateRQPDU](#), [gdcmm::network::AAssociateACPDU](#), [gdcmm::network::PDataTFPDU](#), [gdcmm::network::AAabortPDU](#), [gdcmm::network::AAssociateRJPDU](#), [gdcmm::network::AReleaseRPPDU](#), and [gdcmm::network::AReleaseRQPDU](#).

27.29.3.2 `virtual void gdcmm::network::BasePDU::Print ( std::ostream & os ) const [pure virtual]`

Implemented in [gdcmm::network::AAssociateRQPDU](#), [gdcmm::network::AAssociateACPDU](#), [gdcmm::network::PDataTFPDU](#), [gdcmm::network::AAabortPDU](#), [gdcmm::network::AReleaseRPPDU](#), [gdcmm::network::AReleaseRQPDU](#), and [gdcmm::network::AAssociateRJPDU](#).

27.29.3.3 `virtual std::istream& gdcmm::network::BasePDU::Read ( std::istream & is ) [pure virtual]`

Implemented in [gdcmm::network::AAssociateACPDU](#), [gdcmm::network::AAssociateRQPDU](#), [gdcmm::network::PDataTFPDU](#), [gdcmm::network::AAssociateRJPDU](#), [gdcmm::network::AReleaseRPPDU](#), [gdcmm::network::AReleaseRQPDU](#), and [gdcmm::network::AAabortPDU](#).

27.29.3.4 `virtual size_t gdcmm::network::BasePDU::Size ( ) const [pure virtual]`

Implemented in [gdcmm::network::AAssociateACPDU](#), [gdcmm::network::AAssociateRQPDU](#), [gdcmm::network::PDataTFPDU](#), [gdcmm::network::AAabortPDU](#), [gdcmm::network::AAssociateRJPDU](#), [gdcmm::network::AReleaseRPPDU](#), and [gdcmm::network::AReleaseRQPDU](#).

27.29.3.5 `virtual const std::ostream& gdcmm::network::BasePDU::Write ( std::ostream & os ) const [pure virtual]`

Implemented in [gdcmm::network::AAssociateACPDU](#), [gdcmm::network::AAssociateRQPDU](#), [gdcmm::network::PDataTFPDU](#), [gdcmm::network::AAssociateRJPDU](#), [gdcmm::network::AReleaseRPPDU](#), [gdcmm::network::AReleaseRQPDU](#), and [gdcmm::network::AAabortPDU](#).

The documentation for this class was generated from the following file:


- [gdcmmBasePDU.h](#)

## 27.30 gdcmm::BaseRootQuery Class Reference


[BaseRootQuery](#) contains: a baseclass which will produce a dataset for c-find and c-move with patient/study root.

```
#include <gdcmmBaseRootQuery.h>
```

Inheritance diagram for gdcmm::BaseRootQuery:

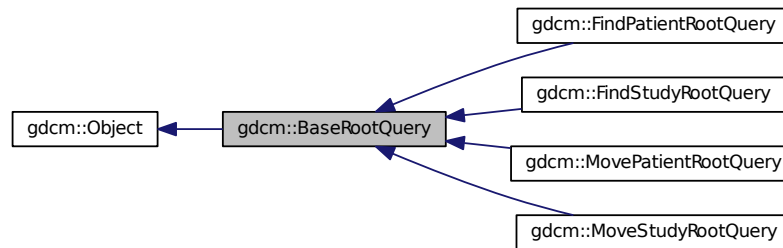

Collaboration diagram for gdcmm::BaseRootQuery:

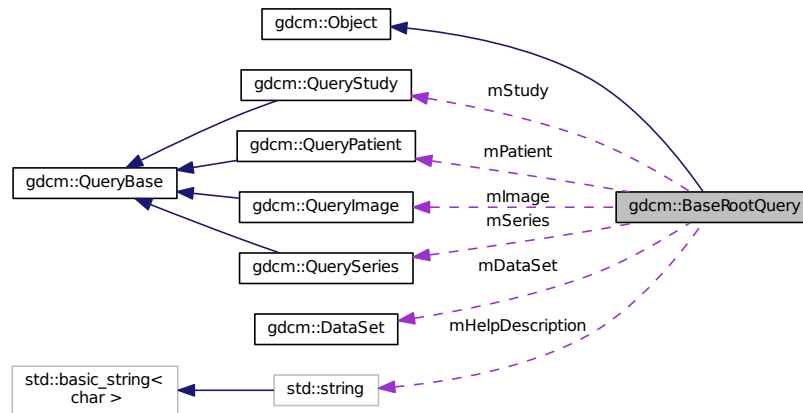

## Public Member Functions

- virtual `~BaseRootQuery ()`
- void `AddQueryDataSet (const DataSet &ds)`
- virtual `UIDs::TSName GetAbstractSyntaxUID () const =0`
- `DataSet const & GetQueryDataSet () const`  
*Set/Get the internal representation of the query as a DataSet.*
- `DataSet & GetQueryDataSet ()`
- `EQueryLevel GetQueryLevelFromQueryRoot (ERootType roottype)`
- virtual `std::vector< Tag > GetTagListByLevel (const EQueryLevel &inQueryLevel)=0`
- virtual void `InitializeDataSet (const EQueryLevel &inQueryLevel)=0`
- void `Print (std::ostream &os) const`
- void `SetSearchParameter (const Tag &inTag, const std::string &inValue)`
- void `SetSearchParameter (const std::string &inKeyword, const std::string &inValue)`
- virtual bool `ValidateQuery (bool inStrict=true) const =0`

- const std::ostream & [WriteHelpFile](#) (std::ostream &os)
- bool [WriteQuery](#) (const std::string &inFileName)

### Static Public Member Functions

- static [QueryBase](#) \* [Construct](#) ([ERootType](#) inRootType, [EQueryLevel](#) qllevel)
- static int [GetQueryLevelFromString](#) (const char \*str)
- static const char \* [GetQueryLevelString](#) ([EQueryLevel](#) ql)

### Protected Member Functions

- [BaseRootQuery](#) ()
- void [SetSearchParameter](#) (const [Tag](#) &inTag, const [DictEntry](#) &inDictEntry, const std::string &inValue)

### Protected Attributes

- [DataSet](#) mDataSet
- std::string mHelpDescription
- [QueryImage](#) mImage
- [QueryPatient](#) mPatient
- [ERootType](#) mRootType
- [QuerySeries](#) mSeries
- [QueryStudy](#) mStudy

### Friends

- class [QueryFactory](#)

## 27.30.1 Detailed Description

[BaseRootQuery](#) contains: a baseclass which will produce a dataset for c-find and c-move with patient/study root.

This class contains the functionality used in patient c-find and c-move queries. PatientRootQuery and StudyRootQuery derive from this class.

Namely: 1) list all tags associated with a particular query type 2) produce a query dataset via tag association

Eventually, it can be used to validate a particular dataset type.

The dataset held by this object (or, really, one of its derivatives) should be passed to a c-find or c-move query.

## 27.30.2 Constructor & Destructor Documentation

27.30.2.1 `gdcm::BaseRootQuery::BaseRootQuery ( )` [protected]

27.30.2.2 `virtual gdcm::BaseRootQuery::~~BaseRootQuery ( )` [virtual]

## 27.30.3 Member Function Documentation

27.30.3.1 void gdcm::BaseRootQuery::AddQueryDataSet ( const DataSet & ds )

27.30.3.2 static QueryBase\* gdcm::BaseRootQuery::Construct ( ERootType inRootType, EQueryLevel qllevel )  
[static]

27.30.3.3 virtual UIDs::TSName gdcm::BaseRootQuery::GetAbstractSyntaxUID ( ) const [pure virtual]

Implemented in [gdcm::FindStudyRootQuery](#), [gdcm::MovePatientRootQuery](#), [gdcm::MoveStudyRootQuery](#), and [gdcm::FindPatientRootQuery](#).

27.30.3.4 DataSet const& gdcm::BaseRootQuery::GetQueryDataSet ( ) const

Set/Get the internal representation of the query as a [DataSet](#).

27.30.3.5 DataSet& gdcm::BaseRootQuery::GetQueryDataSet ( )

27.30.3.6 EQueryLevel gdcm::BaseRootQuery::GetQueryLevelFromQueryRoot ( ERootType roottype )

27.30.3.7 static int gdcm::BaseRootQuery::GetQueryLevelFromString ( const char \* str ) [static]

27.30.3.8 static const char\* gdcm::BaseRootQuery::GetQueryLevelString ( EQueryLevel ql ) [static]

27.30.3.9 virtual std::vector<Tag> gdcm::BaseRootQuery::GetTagListByLevel ( const EQueryLevel & inQueryLevel ) [pure virtual]

this function will return all tags at a given query level, so that they maybe selected for searching. The boolean forFind is true if the query is a find query, or false for a move query.

Implemented in [gdcm::FindPatientRootQuery](#), [gdcm::FindStudyRootQuery](#), [gdcm::MovePatientRootQuery](#), and [gdcm::MoveStudyRootQuery](#).

27.30.3.10 virtual void gdcm::BaseRootQuery::InitializeDataSet ( const EQueryLevel & inQueryLevel ) [pure virtual]

this function sets tag 8,52 to the appropriate value based on query level also fills in the right unique tags, as per the standard's requirements should allow for connection with dcmTk

Implemented in [gdcm::FindPatientRootQuery](#), [gdcm::FindStudyRootQuery](#), [gdcm::MovePatientRootQuery](#), and [gdcm::MoveStudyRootQuery](#).

27.30.3.11 void gdcm::BaseRootQuery::Print ( std::ostream & os ) const [virtual]

Reimplemented from [gdcm::Object](#).

27.30.3.12 void gdcm::BaseRootQuery::SetSearchParameter ( const Tag & inTag, const DictEntry & inDictEntry, const std::string & inValue ) [protected]

27.30.3.13 void gdcm::BaseRootQuery::SetSearchParameter ( const Tag & inTag, const std::string & inValue )

27.30.3.14 void gdcm::BaseRootQuery::SetSearchParameter ( const std::string & inKeyword, const std::string & inValue )

27.30.3.15 `virtual bool gdcm::BaseRootQuery::ValidateQuery ( bool inStrict =true ) const` [pure virtual]

have to be able to ensure that 0x8,0x52 is set (which will be true if InitializeDataSet is called...) that the level is appropriate (ie, not setting PATIENT for a study query that the tags in the query match the right level (either required, unique, optional) by default, this function checks to see if the query is for finding, which is more permissive than for moving. For moving, only the unique tags are allowed. 10 Jan 2011: adding in the 'strict' mode. according to the standard (at least, how I've read it), only tags for a particular level should be allowed in a particular query (ie, just series level tags in a series level query). However, it seems that dcm4chee doesn't share that interpretation. So, if 'inStrict' is false, then tags from the current level and all higher levels are now considered valid. So, if you're doing a non-strict series-level query, tags from the patient and study level can be passed along as well.

Implemented in [gdcm::FindStudyRootQuery](#), [gdcm::MovePatientRootQuery](#), [gdcm::MoveStudyRootQuery](#), and [gdcm::FindPatientRootQuery](#).

27.30.3.16 `const std::ostream& gdcm::BaseRootQuery::WriteHelpFile ( std::ostream & os )`

27.30.3.17 `bool gdcm::BaseRootQuery::WriteQuery ( const std::string & inFileName )`

## 27.30.4 Friends And Related Function Documentation

27.30.4.1 `friend class QueryFactory` [friend]

## 27.30.5 Member Data Documentation

27.30.5.1 `DataSet gdcm::BaseRootQuery::mDataSet` [protected]

27.30.5.2 `std::string gdcm::BaseRootQuery::mHelpDescription` [protected]

27.30.5.3 `QueryImage gdcm::BaseRootQuery::mImage` [protected]

27.30.5.4 `QueryPatient gdcm::BaseRootQuery::mPatient` [protected]

27.30.5.5 `ERootType gdcm::BaseRootQuery::mRootType` [protected]

27.30.5.6 `QuerySeries gdcm::BaseRootQuery::mSeries` [protected]

27.30.5.7 `QueryStudy gdcm::BaseRootQuery::mStudy` [protected]

The documentation for this class was generated from the following file:

- [gdcmBaseRootQuery.h](#)


## 27.31 gdcm::SegmentHelper::BasicCodedEntry Struct Reference

This structure defines a basic coded entry with all of its attributes.

```
#include <gdcmSegmentHelper.h>
```



Collaboration diagram for gdcm::SegmentHelper::BasicCodedEntry:

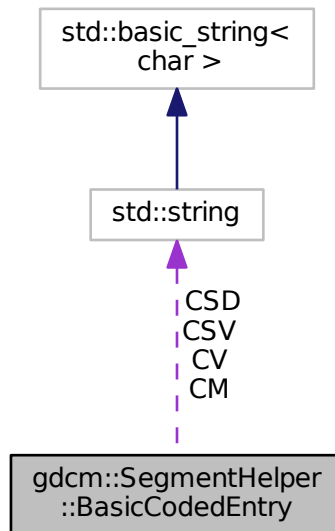

## Public Member Functions

- [BasicCodedEntry](#) ()  
*Constructor.*
- [BasicCodedEntry](#) (const char \*a\_CV, const char \*a\_CSD, const char \*a\_CM)  
*constructor which defines type 1 attributes.*
- [BasicCodedEntry](#) (const char \*a\_CV, const char \*a\_CSD, const char \*a\_CSV, const char \*a\_CM)  
*constructor which defines attributes.*
- bool [IsEmpty](#) (const bool checkOptionalAttributes=false) const  
*Check if each attributes of the basic coded entry is defined.*

## Public Attributes

- std::string [CM](#)  
*Coding Scheme [Version](#) attribute.*
- std::string [CSD](#)  
*Code [Value](#) attribute.*
- std::string [CSV](#)  
*Coding Scheme Designator attribute.*
- std::string [CV](#)

### 27.31.1 Detailed Description

This structure defines a basic coded entry with all of its attributes.

See also

PS 3.3 section 8.8.

### 27.31.2 Constructor & Destructor Documentation

27.31.2.1 `gdcm::SegmentHelper::BasicCodedEntry::BasicCodedEntry ( ) [inline]`

Constructor.

27.31.2.2 `gdcm::SegmentHelper::BasicCodedEntry::BasicCodedEntry ( const char * a_CV, const char * a_CSD, const char * a_CM ) [inline]`

constructor which defines type 1 attributes.

27.31.2.3 `gdcm::SegmentHelper::BasicCodedEntry::BasicCodedEntry ( const char * a_CV, const char * a_CSD, const char * a_CSV, const char * a_CM ) [inline]`

constructor which defines attributes.

### 27.31.3 Member Function Documentation

27.31.3.1 `bool gdcm::SegmentHelper::BasicCodedEntry::IsEmpty ( const bool checkOptionalAttributes = false ) const`

Check if each attributes of the basic coded entry is defined.

Parameters

<i>checkOptionalAttributes</i>	Check also type 1C attributes.
--------------------------------	--------------------------------

### 27.31.4 Member Data Documentation

27.31.4.1 `std::string gdcm::SegmentHelper::BasicCodedEntry::CM`

Coding Scheme [Version](#) attribute.

27.31.4.2 `std::string gdcm::SegmentHelper::BasicCodedEntry::CSD`

Code [Value](#) attribute.

27.31.4.3 `std::string gdcm::SegmentHelper::BasicCodedEntry::CSV`

Coding Scheme Designator attribute.

27.31.4.4 `std::string gdcm::SegmentHelper::BasicCodedEntry::CV`

The documentation for this struct was generated from the following file:


- [gdcmSegmentHelper.h](#)

## 27.32 gdcm::BasicOffsetTable Class Reference


Class to represent a [BasicOffsetTable](#).

```
#include <gdcmBasicOffsetTable.h>
```

Inheritance diagram for gdcm::BasicOffsetTable:

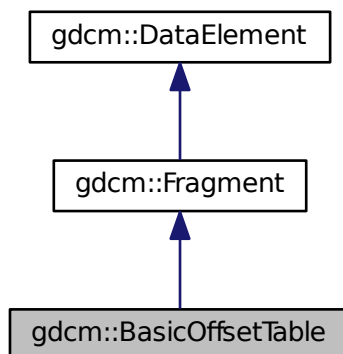

Collaboration diagram for `gdcm::BasicOffsetTable`:

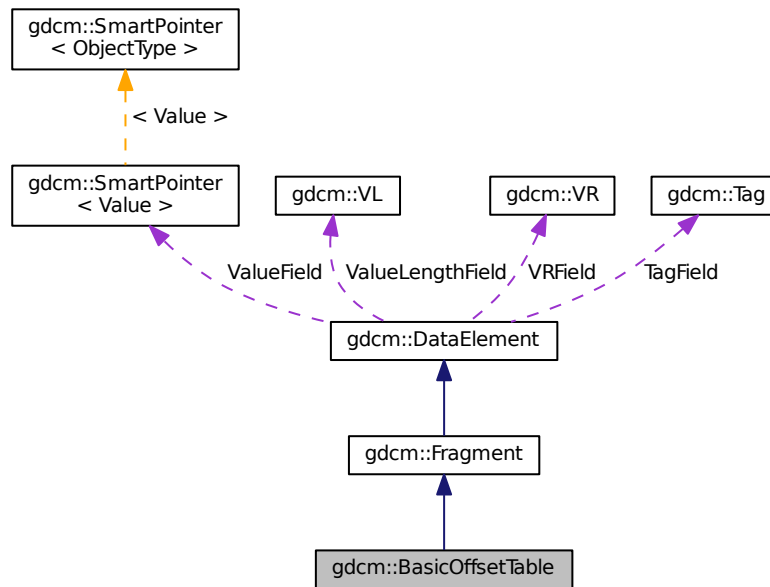

## Public Member Functions

- [BasicOffsetTable](#) ()
- `template<typename TSwap >`  
`std::istream & Read (std::istream &is)`

## Friends

- `std::ostream & operator<< (std::ostream &os, const BasicOffsetTable &val)`

## Additional Inherited Members

### 27.32.1 Detailed Description

Class to represent a [BasicOffsetTable](#).

### 27.32.2 Constructor & Destructor Documentation

#### 27.32.2.1 `gdcm::BasicOffsetTable::BasicOffsetTable ( )` `[inline]`

### 27.32.3 Member Function Documentation

27.32.3.1 `template<typename TSwap> std::istream& gdcm::BasicOffsetTable::Read ( std::istream & is ) [inline]`

References `gdcmDebugMacro`.

## 27.32.4 Friends And Related Function Documentation

27.32.4.1 `std::ostream& operator<< ( std::ostream & os, const BasicOffsetTable & val ) [friend]`

The documentation for this class was generated from the following file:


- [gdcmBasicOffsetTable.h](#)

## 27.33 gdcm::Bitmap Class Reference


**Bitmap** class A bitmap based image. Used as parent for both `IconImage` and the main Pixel Data `Image` It does not contains any World Space information (IPP, IOP)

```
#include <gdcmBitmap.h>
```

Inheritance diagram for `gdcm::Bitmap`:

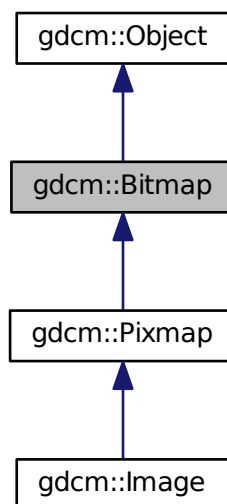

Collaboration diagram for `gdcm::Bitmap`:

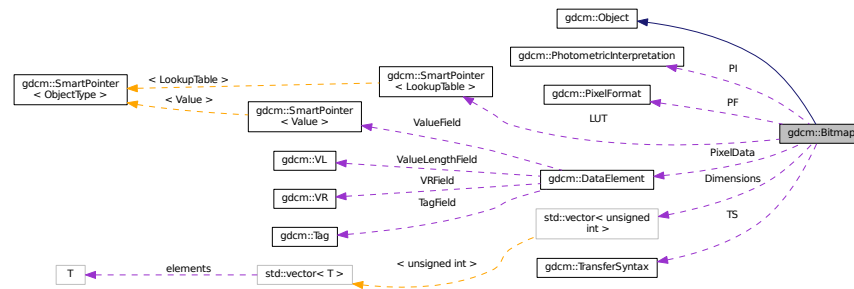

## Public Member Functions

- [Bitmap](#) ()
- [~Bitmap](#) ()
- virtual bool [AreOverlaysInPixelData](#) () const
- void [Clear](#) ()
- bool [GetBuffer](#) (char \*buffer) const  
*Acces the raw data.*
- unsigned long [GetBufferLength](#) () const
- unsigned int [GetColumns](#) () const
- const [DataElement](#) & [GetDataElement](#) () const
- [DataElement](#) & [GetDataElement](#) ()
- unsigned int [GetDimension](#) (unsigned int idx) const
- const unsigned int \* [GetDimensions](#) () const  
*Return the dimension of the pixel data, first dimension (x), then 2nd (y), then 3rd (z)...*
- const [LookupTable](#) & [GetLUT](#) () const
- [LookupTable](#) & [GetLUT](#) ()
- bool [GetNeedByteSwap](#) () const
- unsigned int [GetNumberOfDimensions](#) () const  
*Return the number of dimension of the pixel data bytes; for example 2 for a 2D matrices of values.*
- const [PhotometricInterpretation](#) & [GetPhotometricInterpretation](#) () const  
*return the photometric interpretation*
- const [PixelFormat](#) & [GetPixelFormat](#) () const  
*Get/Set PixelFormat.*
- [PixelFormat](#) & [GetPixelFormat](#) ()
- unsigned int [GetPlanarConfiguration](#) () const  
*return the planar configuration*
- unsigned int [GetRows](#) () const
- const [TransferSyntax](#) & [GetTransferSyntax](#) () const
- bool [IsEmpty](#) () const
- bool [IsLossy](#) () const  
*Return whether or not the image was compressed using a lossy compressor or not.*
- bool [IsTransferSyntaxCompatible](#) ([TransferSyntax](#) const &ts) const
- void [Print](#) (std::ostream &) const
- void [SetColumns](#) (unsigned int col)

- void [SetDataElement](#) ([DataElement](#) const &de)
- void [SetDimension](#) (unsigned int idx, unsigned int dim)
- void [SetDimensions](#) (const unsigned int dims[3])
- void [SetLossyFlag](#) (bool f)  
Specifically set that the image was compressed using a lossy compression mechanism.
- void [SetLUT](#) ([LookupTable](#) const &lut)  
Set/Get LUT.
- void [SetNeedByteSwap](#) (bool b)
- void [SetNumberOfDimensions](#) (unsigned int dim)
- void [SetPhotometricInterpretation](#) ([PhotometricInterpretation](#) const &pi)
- void [SetPixelFormat](#) ([PixelFormat](#) const &pf)
- void [SetPlanarConfiguration](#) (unsigned int pc)
- void [SetRows](#) (unsigned int rows)
- void [SetTransferSyntax](#) ([TransferSyntax](#) const &ts)  
Transfer syntax.

## Protected Types

- typedef [SmartPointer](#)< [LookupTable](#) > [LUTPtr](#)

## Protected Member Functions

- bool [ComputeLossyFlag](#) ()
- bool [GetBuffer2](#) (std::ostream &os) const
- bool [TryJPEG2000Codec](#) (char \*buffer, bool &lossyflag) const
- bool [TryJPEG2000Codec2](#) (std::ostream &os) const
- bool [TryJPEGCodec](#) (char \*buffer, bool &lossyflag) const
- bool [TryJPEGCodec2](#) (std::ostream &os) const
- bool [TryJPEGLSCodec](#) (char \*buffer, bool &lossyflag) const
- bool [TryKAKADUCodec](#) (char \*buffer, bool &lossyflag) const
- bool [TryPVRGCodec](#) (char \*buffer, bool &lossyflag) const
- bool [TryRAWCodec](#) (char \*buffer, bool &lossyflag) const
- bool [TryRLECodec](#) (char \*buffer, bool &lossyflag) const

## Protected Attributes

- std::vector< unsigned int > [Dimensions](#)
- bool [LossyFlag](#)
- [LUTPtr](#) [LUT](#)
- bool [NeedByteSwap](#)
- unsigned int [NumberOfDimensions](#)
- [PixelFormat](#) [PF](#)
- [PhotometricInterpretation](#) [PI](#)
- [DataElement](#) [PixelData](#)
- unsigned int [PlanarConfiguration](#)
- [TransferSyntax](#) [TS](#)

## Friends

- class [ImageChangeTransferSyntax](#)
- class [PixmapReader](#)

### 27.33.1 Detailed Description

[Bitmap](#) class A bitmap based image. Used as parent for both [IconImage](#) and the main Pixel Data [Image](#) It does not contains any World Space information (IPP, IOP)

Examples:

[ExtractIconFromFile.cxx](#).

### 27.33.2 Member Typedef Documentation

27.33.2.1 `typedef SmartPointer<LookupTable> gdcm::Bitmap::LUTPtr` `[protected]`

### 27.33.3 Constructor & Destructor Documentation

27.33.3.1 `gdcm::Bitmap::Bitmap ( )`

27.33.3.2 `gdcm::Bitmap::~~Bitmap ( )`

### 27.33.4 Member Function Documentation

27.33.4.1 `virtual bool gdcm::Bitmap::AreOverlaysInPixelData ( ) const` `[inline],[virtual]`

Reimplemented in [gdcm::Pixmap](#).

27.33.4.2 `void gdcm::Bitmap::Clear ( )`

27.33.4.3 `bool gdcm::Bitmap::ComputeLossyFlag ( )` `[protected]`

27.33.4.4 `bool gdcm::Bitmap::GetBuffer ( char * buffer ) const`

Acces the raw data.

Examples:

[ConvertToQImage.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

27.33.4.5 `bool gdcm::Bitmap::GetBuffer2 ( std::ostream & os ) const` `[protected]`

27.33.4.6 `unsigned long gdcm::Bitmap::GetBufferLength ( ) const`

Return the length of the image after decompression WARNING for palette color: It will NOT take into account the Palette Color thus you need to multiply this length by 3 if computing the size of equivalent RGB image



Examples:

[ConvertToQImage.cxx](#), [GenFakelImage.cxx](#), [GetSubSequenceData.cxx](#), [PatchFile.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

27.33.4.7 `unsigned int gdcm::Bitmap::GetColumns ( ) const [inline]`

27.33.4.8 `const DataElement& gdcm::Bitmap::GetDataElement ( ) const [inline]`

Examples:

[ExtractIconFromFile.cxx](#).

27.33.4.9 `DataElement& gdcm::Bitmap::GetDataElement ( ) [inline]`

27.33.4.10 `unsigned int gdcm::Bitmap::GetDimension ( unsigned int idx ) const`

27.33.4.11 `const unsigned int* gdcm::Bitmap::GetDimensions ( ) const`

Return the dimension of the pixel data, first dimension (x), then 2nd (y), then 3rd (z)...

Examples:

[ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), [FixJAIBugJPEGLS.cxx](#), [HelloVizWorld.cxx](#), and [threadgdcm.cxx](#).

27.33.4.12 `const LookupTable& gdcm::Bitmap::GetLUT ( ) const [inline]`

Examples:

[ExtractIconFromFile.cxx](#).

27.33.4.13 `LookupTable& gdcm::Bitmap::GetLUT ( ) [inline]`

27.33.4.14 `bool gdcm::Bitmap::GetNeedByteSwap ( ) const [inline]`

27.33.4.15 `unsigned int gdcm::Bitmap::GetNumberOfDimensions ( ) const`

Return the number of dimension of the pixel data bytes; for example 2 for a 2D matrices of values.

Examples:

[HelloVizWorld.cxx](#), and [threadgdcm.cxx](#).

27.33.4.16 `const PhotometricInterpretation& gdcm::Bitmap::GetPhotometricInterpretation ( ) const`

return the photometric interpretation

Examples:

[ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), and [HelloVizWorld.cxx](#).

27.33.4.17 `const PixelFormat& gdcm::Bitmap::GetPixelFormat ( ) const` `[inline]`

Get/Set [PixelFormat](#).

Examples:

[ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenFakelImage.cxx](#), [GetJPEGSamplePrecision.cxx](#), and [threadgdcm.cxx](#).

27.33.4.18 `PixelFormat& gdcm::Bitmap::GetPixelFormat ( )` `[inline]`

27.33.4.19 `unsigned int gdcm::Bitmap::GetPlanarConfiguration ( ) const`

return the planar configuration

27.33.4.20 `unsigned int gdcm::Bitmap::GetRows ( ) const` `[inline]`

27.33.4.21 `const TransferSyntax& gdcm::Bitmap::GetTransferSyntax ( ) const` `[inline]`

Examples:

[ExtractIconFromFile.cxx](#).

27.33.4.22 `bool gdcm::Bitmap::IsEmpty ( ) const` `[inline]`

27.33.4.23 `bool gdcm::Bitmap::IsLossy ( ) const`

Return whether or not the image was compressed using a lossy compressor or not.

27.33.4.24 `bool gdcm::Bitmap::IsTransferSyntaxCompatible ( TransferSyntax const & ts ) const`

27.33.4.25 `void gdcm::Bitmap::Print ( std::ostream & ) const` `[virtual]`

Reimplemented from [gdcm::Object](#).

Reimplemented in [gdcm::Image](#), and [gdcm::Pixmap](#).

Examples:

[ExtractIconFromFile.cxx](#).

27.33.4.26 `void gdcm::Bitmap::SetColumns ( unsigned int col )` `[inline]`

27.33.4.27 `void gdcm::Bitmap::SetDataElement ( DataElement const & de )` `[inline]`

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [GenFakelImage.cxx](#), [GetSubSequenceData.cxx](#), and [iU22tomultisc.cxx](#).

27.33.4.28 void gdcm::Bitmap::SetDimension ( unsigned int *idx*, unsigned int *dim* )

Examples:

[csa2img.cxx](#), [GenFakelImage.cxx](#), [GetSubSequenceData.cxx](#), and [iU22tomultisc.cxx](#).

27.33.4.29 void gdcm::Bitmap::SetDimensions ( const unsigned int *dims*[3] )

Examples:

[CreateARGBImage.cxx](#), and [CreateCMYKImage.cxx](#).

27.33.4.30 void gdcm::Bitmap::SetLossyFlag ( bool *f* ) [inline]

Specifically set that the image was compressed using a lossy compression mechanism.

27.33.4.31 void gdcm::Bitmap::SetLUT ( LookupTable const & *lut* ) [inline]

Set/Get LUT.

27.33.4.32 void gdcm::Bitmap::SetNeedByteSwap ( bool *b* ) [inline]

27.33.4.33 void gdcm::Bitmap::SetNumberOfDimensions ( unsigned int *dim* )

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [FileChangeTS.cs](#), [FileChangeTSLossy.cs](#), [GenFakelImage.cxx](#), [GetSubSequenceData.cxx](#), and [iU22tomultisc.cxx](#).

27.33.4.34 void gdcm::Bitmap::SetPhotometricInterpretation ( PhotometricInterpretation const & *pi* )

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [GenFakelImage.cxx](#), [GetSubSequenceData.cxx](#), and [iU22tomultisc.cxx](#).

27.33.4.35 void gdcm::Bitmap::SetPixelFormat ( PixelFormat const & *pf* ) [inline]

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), and [iU22tomultisc.cxx](#).

References gdcm::PixelFormat::Validate().

27.33.4.36 void gdcm::Bitmap::SetPlanarConfiguration ( unsigned int *pc* )

Warning

you need to call SetPixelFormat first (before SetPlanarConfiguration) for consistency checking

27.33.4.37 void gdcmm::Bitmap::SetRows ( unsigned int *rows* ) [inline]

27.33.4.38 void gdcmm::Bitmap::SetTransferSyntax ( TransferSyntax const & *ts* ) [inline]

Transfer syntax.

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), and [MergeTwoFiles.cxx](#).

27.33.4.39 bool gdcmm::Bitmap::TryJPEG2000Codec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.40 bool gdcmm::Bitmap::TryJPEG2000Codec2 ( std::ostream & *os* ) const [protected]

27.33.4.41 bool gdcmm::Bitmap::TryJPEGCodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.42 bool gdcmm::Bitmap::TryJPEGCodec2 ( std::ostream & *os* ) const [protected]

27.33.4.43 bool gdcmm::Bitmap::TryJPEGLSCodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.44 bool gdcmm::Bitmap::TryKAKADUCodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.45 bool gdcmm::Bitmap::TryPVRGCodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.46 bool gdcmm::Bitmap::TryRAWCodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

27.33.4.47 bool gdcmm::Bitmap::TryRLECodec ( char \* *buffer*, bool & *lossyflag* ) const [protected]

## 27.33.5 Friends And Related Function Documentation

27.33.5.1 friend class ImageChangeTransferSyntax [friend]

27.33.5.2 friend class PixmapReader [friend]

## 27.33.6 Member Data Documentation

27.33.6.1 std::vector<unsigned int> gdcmm::Bitmap::Dimensions [protected]

27.33.6.2 bool gdcmm::Bitmap::LossyFlag [protected]

27.33.6.3 LUTPtr gdcmm::Bitmap::LUT [protected]

27.33.6.4 bool gdcmm::Bitmap::NeedByteSwap [protected]

27.33.6.5 unsigned int gdcmm::Bitmap::NumberOfDimensions [protected]

27.33.6.6 PixelFormat gdcmm::Bitmap::PF [protected]

27.33.6.7 PhotometricInterpretation gdcmm::Bitmap::PI [protected]

27.33.6.8 DataElement gdcmm::Bitmap::PixelData [protected]

27.33.6.9 unsigned int gdcm::Bitmap::PlanarConfiguration [protected]

27.33.6.10 TransferSyntax gdcm::Bitmap::TS [protected]

The documentation for this class was generated from the following file:


- [gdcmBitmap.h](#)

## 27.34 gdcm::BitmapToBitmapFilter Class Reference


[BitmapToBitmapFilter](#) class Super class for all filter taking an image and producing an output image.

```
#include <gdcmBitmapToBitmapFilter.h>
```

Inheritance diagram for gdcm::BitmapToBitmapFilter:

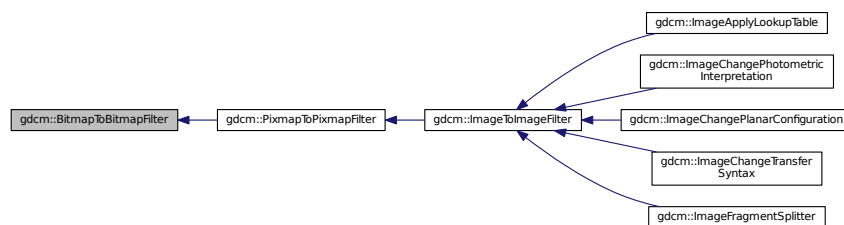

Collaboration diagram for gdcm::BitmapToBitmapFilter:

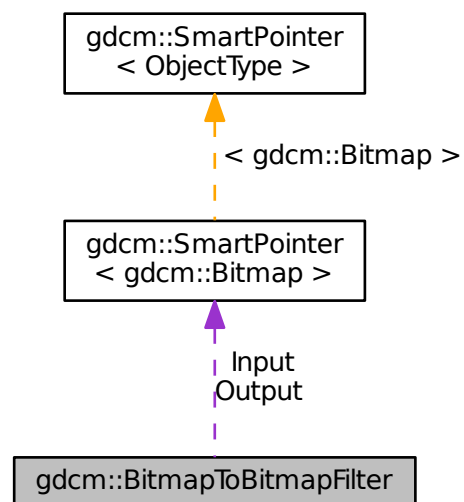

## Public Member Functions

- [BitmapToBitmapFilter](#) ()
- [~BitmapToBitmapFilter](#) ()
- const [Bitmap](#) & [GetOutput](#) () const  
*Get Output image.*
- const [Bitmap](#) & [GetOutputAsBitmap](#) () const
- void [SetInput](#) (const [Bitmap](#) &image)  
*Set input image.*

## Protected Attributes

- [SmartPointer](#)< [Bitmap](#) > [Input](#)
- [SmartPointer](#)< [Bitmap](#) > [Output](#)

### 27.34.1 Detailed Description

[BitmapToBitmapFilter](#) class Super class for all filter taking an image and producing an output image.

### 27.34.2 Constructor & Destructor Documentation

27.34.2.1 `gdcm::BitmapToBitmapFilter::BitmapToBitmapFilter ( )`

27.34.2.2 `gdcm::BitmapToBitmapFilter::~~BitmapToBitmapFilter ( )` `[inline]`

### 27.34.3 Member Function Documentation

27.34.3.1 `const Bitmap& gdcm::BitmapToBitmapFilter::GetOutput ( ) const` `[inline]`

Get Output image.

27.34.3.2 `const Bitmap& gdcm::BitmapToBitmapFilter::GetOutputAsBitmap ( ) const`

27.34.3.3 `void gdcm::BitmapToBitmapFilter::SetInput ( const Bitmap & image )`

Set input image.

Examples:

[CompressImage.cxx](#).

### 27.34.4 Member Data Documentation

27.34.4.1 `SmartPointer<Bitmap> gdcm::BitmapToBitmapFilter::Input` `[protected]`

27.34.4.2 `SmartPointer<Bitmap> gdcm::BitmapToBitmapFilter::Output` `[protected]`

The documentation for this class was generated from the following file:


- [gdcmBitmapToBitmapFilter.h](#)

## 27.35 gdcm::BoxRegion Class Reference


Class for manipulation box region This is a very simple implementation of the [Region](#) class. It only support 3D box type region. It assumes the 3D Box does not have a tilt Origin is as (0,0,0)

```
#include <gdcmBoxRegion.h>
```

Inheritance diagram for gdcm::BoxRegion:

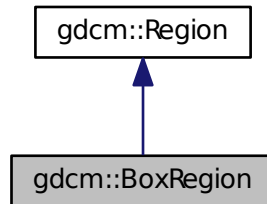

Collaboration diagram for gdcm::BoxRegion:

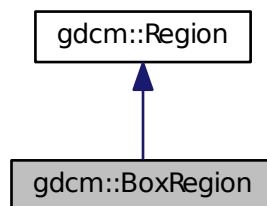

### Public Member Functions

- [BoxRegion](#) ()
- [BoxRegion](#) (const [BoxRegion](#) &)  
*copy/cstor and al.*
- [~BoxRegion](#) ()
- [size\\_t Area](#) () const  
*compute the area*
- [Region \\* Clone](#) () const
- [BoxRegion ComputeBoundingBox](#) ()  
*Return the Axis-Aligned minimum bounding box for all regions.*

- bool [Empty](#) () const  
*return whether this domain is empty:*
- unsigned int [GetXMax](#) () const
- unsigned int [GetXMin](#) () const  
*Get domain.*
- unsigned int [GetYMax](#) () const
- unsigned int [GetYMin](#) () const
- unsigned int [GetZMax](#) () const
- unsigned int [GetZMin](#) () const
- bool [IsValid](#) () const  
*return whether this is valid domain*
- void [operator=](#) (const [BoxRegion](#) &)
- void [Print](#) (std::ostream &os=std::cout) const  
*Print.*
- void [SetDomain](#) (unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax)  
*Set domain.*

## Static Public Member Functions

- static [BoxRegion BoundingBox](#) ([BoxRegion](#) const &b1, [BoxRegion](#) const &b2)  
*Helper class to compute the bounding box of two [BoxRegion](#).*

## 27.35.1 Detailed Description

Class for manipulation box region This is a very simple implementation of the [Region](#) class. It only support 3D box type region. It assumes the 3D Box does not have a tilt Origin is as (0,0,0)

## 27.35.2 Constructor & Destructor Documentation

27.35.2.1 `gdcm::BoxRegion::BoxRegion ( )`

27.35.2.2 `gdcm::BoxRegion::~~BoxRegion ( )`

27.35.2.3 `gdcm::BoxRegion::BoxRegion ( const BoxRegion & )`

copy/cstor and al.

## 27.35.3 Member Function Documentation

27.35.3.1 `size_t gdcm::BoxRegion::Area ( ) const` `[virtual]`

compute the area

Implements [gdcm::Region](#).



27.35.3.2 **static BoxRegion** gdcm::BoxRegion::BoundingBox ( **BoxRegion** const & *b1*, **BoxRegion** const & *b2* )  
[static]

Helper class to compute the bounding box of two [BoxRegion](#).

27.35.3.3 **Region\*** gdcm::BoxRegion::Clone ( ) const [virtual]

Implements [gdcm::Region](#).

27.35.3.4 **BoxRegion** gdcm::BoxRegion::ComputeBoundingBox ( ) [virtual]

Return the Axis-Aligned minimum bounding box for all regions.

Implements [gdcm::Region](#).

27.35.3.5 **bool** gdcm::BoxRegion::Empty ( ) const [virtual]

return whether this domain is empty:

Implements [gdcm::Region](#).

27.35.3.6 **unsigned int** gdcm::BoxRegion::GetXMax ( ) const

27.35.3.7 **unsigned int** gdcm::BoxRegion::GetXMin ( ) const

Get domain.

27.35.3.8 **unsigned int** gdcm::BoxRegion::GetYMax ( ) const

27.35.3.9 **unsigned int** gdcm::BoxRegion::GetYMin ( ) const

27.35.3.10 **unsigned int** gdcm::BoxRegion::GetZMax ( ) const

27.35.3.11 **unsigned int** gdcm::BoxRegion::GetZMin ( ) const

27.35.3.12 **bool** gdcm::BoxRegion::IsValid ( ) const [virtual]

return whether this is valid domain

Implements [gdcm::Region](#).

27.35.3.13 **void** gdcm::BoxRegion::operator= ( **const BoxRegion** & )

27.35.3.14 **void** gdcm::BoxRegion::Print ( **std::ostream** & *os* = **std::cout** ) const [virtual]

Print.

Reimplemented from [gdcm::Region](#).

27.35.3.15 void gdcM::BoxRegion::SetDomain ( unsigned int *xmin*, unsigned int *xmax*, unsigned int *ymin*, unsigned int *ymax*, unsigned int *zmin*, unsigned int *zmax* )

Set domain.

The documentation for this class was generated from the following file:

- [gdcMBoxRegion.h](#)

## 27.36 gdcM::ByteBuffer Class Reference

[ByteBuffer](#).

```
#include <gdcMByteBuffer.h>
```

### Public Member Functions

- [ByteBuffer](#) ()
- char \* [Get](#) (int len)
- const char \* [GetStart](#) () const
- void [ShiftEnd](#) (int len)
- void [UpdatePosition](#) ()

### 27.36.1 Detailed Description

[ByteBuffer](#).

Detailed description here

#### Note

looks like a std::streambuf or std::filebuf class with the get and peek pointer

### 27.36.2 Constructor & Destructor Documentation

27.36.2.1 gdcM::ByteBuffer::ByteBuffer ( ) [inline]

### 27.36.3 Member Function Documentation

27.36.3.1 char\* gdcM::ByteBuffer::Get ( int len ) [inline]

27.36.3.2 const char\* gdcM::ByteBuffer::GetStart ( ) const [inline]

27.36.3.3 void gdcM::ByteBuffer::ShiftEnd ( int len ) [inline]

27.36.3.4 void gdcM::ByteBuffer::UpdatePosition ( ) [inline]

The documentation for this class was generated from the following file:

- [gdcMByteBuffer.h](#)

## 27.37 gdcmm::ByteSwap< T > Class Template Reference

[ByteSwap.](#)

```
#include <gdcmmByteSwap.h>
```

### Static Public Member Functions

- static void [Swap](#) (T &p)
- static void [SwapFromSwapCodeIntoSystem](#) (T &p, [SwapCode](#) const &sc)
- static void [SwapRange](#) (T \*p, unsigned int num)
- static void [SwapRangeFromSwapCodeIntoSystem](#) (T \*p, [SwapCode](#) const &sc, std::streamoff num)
- static bool [SystemIsBigEndian](#) ()
- static bool [SystemIsLittleEndian](#) ()

### 27.37.1 Detailed Description

```
template<class T>class gdcmm::ByteSwap< T >
```

[ByteSwap.](#)

Perform machine dependent byte swapping (Little Endian, Big Endian, Bad Little Endian, Bad Big Endian). TODO: bswap\_32 / bswap\_64 ...

Examples:

[TestByteSwap.cxx.](#)

### 27.37.2 Member Function Documentation

27.37.2.1 `template<class T> static void gdcmm::ByteSwap< T >::Swap ( T & p ) [static]`

27.37.2.2 `template<class T> static void gdcmm::ByteSwap< T >::SwapFromSwapCodeIntoSystem ( T & p, SwapCode const & sc ) [static]`

Examples:

[TestByteSwap.cxx.](#)

27.37.2.3 `template<class T> static void gdcmm::ByteSwap< T >::SwapRange ( T * p, unsigned int num ) [static]`

27.37.2.4 `template<class T> static void gdcmm::ByteSwap< T >::SwapRangeFromSwapCodeIntoSystem ( T * p, SwapCode const & sc, std::streamoff num ) [static]`

Examples:

[TestByteSwap.cxx.](#)

27.37.2.5 `template<class T> static bool gdcmm::ByteSwap< T >::SystemIsBigEndian ( ) [static]`

Query the machine Endian-ness.

27.37.2.6 `template<class T> static bool gdcm::ByteSwap<T>::SystemIsLittleEndian ( ) [static]`

The documentation for this class was generated from the following file:

- [gdcmByteSwap.h](#)

## 27.38 gdcm::ByteSwapFilter Class Reference

[ByteSwapFilter](#) In place byte-swapping of a dataset FIXME: FL status ??

```
#include <gdcmByteSwapFilter.h>
```

### Public Member Functions

- [ByteSwapFilter](#) ([DataSet](#) &ds)
- [~ByteSwapFilter](#) ()
- bool [ByteSwap](#) ()
- void [SetByteSwapTag](#) (bool b)

### 27.38.1 Detailed Description

[ByteSwapFilter](#) In place byte-swapping of a dataset FIXME: FL status ??

### 27.38.2 Constructor & Destructor Documentation

27.38.2.1 `gdcm::ByteSwapFilter::ByteSwapFilter ( DataSet & ds ) [inline]`

27.38.2.2 `gdcm::ByteSwapFilter::~~ByteSwapFilter ( )`

### 27.38.3 Member Function Documentation

27.38.3.1 `bool gdcm::ByteSwapFilter::ByteSwap ( )`

27.38.3.2 `void gdcm::ByteSwapFilter::SetByteSwapTag ( bool b ) [inline]`

The documentation for this class was generated from the following file:


- [gdcmByteSwapFilter.h](#)

## 27.39 gdcm::ByteValue Class Reference


Class to represent binary value (array of bytes)

```
#include <gdcmByteValue.h>
```

Inheritance diagram for gdcm::ByteValue:

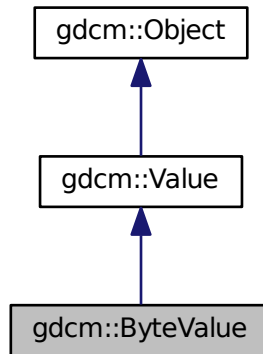

Collaboration diagram for gdcm::ByteValue:

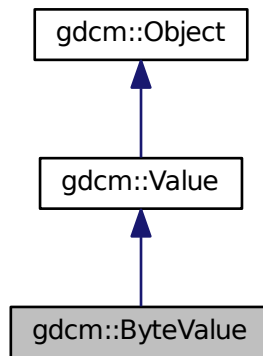

### Public Member Functions

- [ByteValue](#) (const char \*array=0, [VL](#) const &vl=0)
- [ByteValue](#) (std::vector< char > &v)
- [~ByteValue](#) ()
- void [Append](#) ([ByteValue](#) const &bv)
- void [Clear](#) ()
- [VL ComputeLength](#) () const
- void [Fill](#) (char c)

- bool [GetBuffer](#) (char \*buffer, unsigned long length) const
- [VL GetLength](#) () const
- const char \* [GetPointer](#) () const
- bool [IsEmpty](#) () const
- bool [IsPrintable](#) (VL length) const

*Checks whether a 'ByteValue' is printable or not (in order to avoid corrupting the terminal of invocation when printing) / don't think this function is working since it does not handle UNICODE or character set...*

- [operator const std::vector< char > & \(\)](#) const
- [ByteValue & operator=](#) (const [ByteValue](#) &val)
- bool [operator==](#) (const [ByteValue](#) &val) const
- bool [operator==](#) (const [Value](#) &val) const
- void [PrintASCII](#) (std::ostream &os, [VL](#) maxlength) const
- void [PrintASCIIXML](#) (std::ostream &os) const
- void [PrintGroupLength](#) (std::ostream &os)
- void [PrintHex](#) (std::ostream &os, [VL](#) maxlength) const
- void [PrintHexXML](#) (std::ostream &os) const
- void [PrintPNXML](#) (std::ostream &os) const
- template<typename TSwap, typename TType >  
std::istream & [Read](#) (std::istream &is, bool readvalues=true)
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- void [SetLength](#) (VL vl)
- template<typename TSwap, typename TType >  
std::ostream const & [Write](#) (std::ostream &os) const
- template<typename TSwap >  
std::ostream const & [Write](#) (std::ostream &os) const
- bool [WriteBuffer](#) (std::ostream &os) const

## Protected Member Functions

- void [Print](#) (std::ostream &os) const
- void [SetLengthOnly](#) (VL vl)

### 27.39.1 Detailed Description

Class to represent binary value (array of bytes)

Note

Examples:

[DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetSubSequenceData.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), and [rle2img.cxx](#).

### 27.39.2 Constructor & Destructor Documentation

27.39.2.1 `gdcmm::ByteValue::ByteValue ( const char * array = 0, VL const & vl = 0 ) [inline]`

References [gdcmmDebugMacro](#).

27.39.2.2 `gdcm::ByteValue::ByteValue ( std::vector< char > & v ) [inline]`

Warning

casting to uint32\_t

27.39.2.3 `gdcm::ByteValue::~ByteValue ( ) [inline]`

### 27.39.3 Member Function Documentation

27.39.3.1 `void gdcm::ByteValue::Append ( ByteValue const & bv )`

27.39.3.2 `void gdcm::ByteValue::Clear ( ) [inline],[virtual]`

Implements [gdcm::Value](#).

27.39.3.3 `VL gdcm::ByteValue::ComputeLength ( ) const [inline]`

Referenced by `gdcm::Fragment::Write()`.

27.39.3.4 `void gdcm::ByteValue::Fill ( char c ) [inline]`

Examples:

[DuplicatePCDE.cxx](#).

27.39.3.5 `bool gdcm::ByteValue::GetBuffer ( char * buffer, unsigned long length ) const`

Examples:

[FixJAIBugJPEGLS.cxx](#).

27.39.3.6 `VL gdcm::ByteValue::GetLength ( ) const [inline],[virtual]`

Implements [gdcm::Value](#).

Examples:

[DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetSubSequenceData.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcm::operator<<()`, `gdcm::SequenceOfFragments::ReadValue()`, `gdcm::Element< VR::OB, VM::VM1_n >::Set()`, `gdcm::Element< TVR, VM::VM1_n >::Set()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValueNoSwap()`, `gdcm::Element< VR::OB, VM::VM1_n >::SetNoSwap()`, `gdcm::Element< TVR, VM::VM1_n >::SetNoSwap()`, and `gdcm::Fragment::Write()`.

27.39.3.7 `const char* gdcm::ByteValue::GetPointer ( ) const [inline]`

Examples:

[DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [GetSubSequenceData.cxx](#), [MrProtocol.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcm::operator<<()`, `gdcm::SequenceOfFragments::ReadValue()`, `gdcm::Element< VR::OB, VM::VM1_n >::Set()`, `gdcm::Element< TVR, VM::VM1_n >::Set()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::SetByteValue()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetByteValueNoSwap()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetByteValueNoSwap()`, `gdcm::Element< VR::OB, VM::VM1_n >::SetNoSwap()`, and `gdcm::Element< TVR, VM::VM1_n >::SetNoSwap()`.

27.39.3.8 `bool gdcm::ByteValue::IsEmpty ( ) const [inline]`

27.39.3.9 `bool gdcm::ByteValue::IsPrintable ( VL length ) const [inline]`

Checks whether a 'ByteValue' is printable or not (in order to avoid corrupting the terminal of invocation when printing) I don't think this function is working since it does not handle UNICODE or character set...

27.39.3.10 `gdcm::ByteValue::operator const std::vector< char > & ( ) const [inline]`

27.39.3.11 `ByteValue& gdcm::ByteValue::operator= ( const ByteValue & val ) [inline]`

27.39.3.12 `bool gdcm::ByteValue::operator== ( const ByteValue & val ) const [inline]`

27.39.3.13 `bool gdcm::ByteValue::operator== ( const Value & val ) const [inline], [virtual]`

Implements [gdcm::Value](#).

27.39.3.14 `void gdcm::ByteValue::Print ( std::ostream & os ) const [inline], [protected], [virtual]`

Reimplemented from [gdcm::Object](#).

27.39.3.15 `void gdcm::ByteValue::PrintASCII ( std::ostream & os, VL maxlength ) const`

27.39.3.16 `void gdcm::ByteValue::PrintASCIIXML ( std::ostream & os ) const`

27.39.3.17 `void gdcm::ByteValue::PrintGroupLength ( std::ostream & os ) [inline]`

27.39.3.18 `void gdcm::ByteValue::PrintHex ( std::ostream & os, VL maxlength ) const`

27.39.3.19 `void gdcm::ByteValue::PrintHexXML ( std::ostream & os ) const`

27.39.3.20 `void gdcm::ByteValue::PrintPXML ( std::ostream & os ) const`

To Print Values in Native DICOM format



27.39.3.21 `template<typename TSwap, typename TType> std::istream& gdcm::ByteValue::Read ( std::istream & is, bool readvalues = true ) [inline]`

27.39.3.22 `template<typename TSwap> std::istream& gdcm::ByteValue::Read ( std::istream & is ) [inline]`

27.39.3.23 `void gdcm::ByteValue::SetLength ( VL vl ) [inline],[virtual]`

Implements [gdcm::Value](#).

References [gdcmDebugMacro](#), [gdcm::VL::IsOdd\(\)](#), and [gdcm::VL::IsUndefined\(\)](#).

27.39.3.24 `void gdcm::ByteValue::SetLengthOnly ( VL vl ) [inline],[protected],[virtual]`

Reimplemented from [gdcm::Value](#).

27.39.3.25 `template<typename TSwap, typename TType> std::ostream const& gdcm::ByteValue::Write ( std::ostream & os ) const [inline]`

Referenced by [gdcm::Fragment::Write\(\)](#).

27.39.3.26 `template<typename TSwap> std::ostream const& gdcm::ByteValue::Write ( std::ostream & os ) const [inline]`

27.39.3.27 `bool gdcm::ByteValue::WriteBuffer ( std::ostream & os ) const [inline]`


The documentation for this class was generated from the following file:

- [gdcmByteValue.h](#)


## 27.40 gdcm::CAPICryptoFactory Class Reference

```
#include <gdcmCAPICryptoFactory.h>
```

Inheritance diagram for [gdcm::CAPICryptoFactory](#):

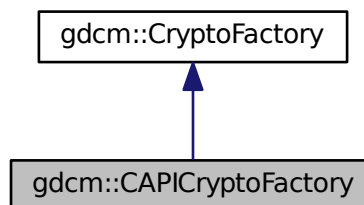

Collaboration diagram for `gdcm::CAPICryptoFactory`:

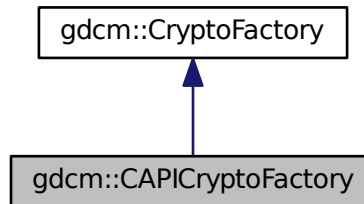

## Public Member Functions

- [CAPICryptoFactory](#) ([CryptoLib](#) id)
- [CryptographicMessageSyntax](#) \* [CreateCMSProvider](#) ()

## Additional Inherited Members

### 27.40.1 Constructor & Destructor Documentation

27.40.1.1 `gdcm::CAPICryptoFactory::CAPICryptoFactory ( CryptoLib id )`

### 27.40.2 Member Function Documentation

27.40.2.1 `CryptographicMessageSyntax* gdcm::CAPICryptoFactory::CreateCMSProvider ( )` [[virtual](#)]

Implements [gdcm::CryptoFactory](#).


The documentation for this class was generated from the following file:

- [gdcmCAPICryptoFactory.h](#)


## 27.41 `gdcm::CAPICryptographicMessageSyntax` Class Reference

```
#include <gdcmCAPICryptographicMessageSyntax.h>
```

Inheritance diagram for gdcM::CAPICryptographicMessageSyntax:

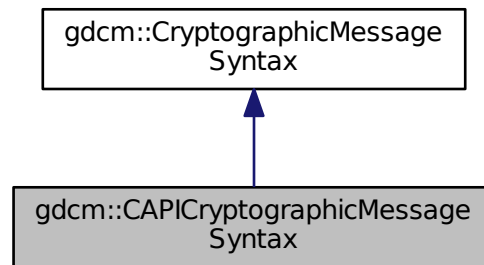

Collaboration diagram for gdcM::CAPICryptographicMessageSyntax:

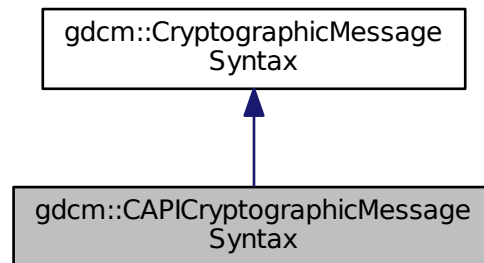

## Public Member Functions

- [CAPICryptographicMessageSyntax](#) ()
- [~CAPICryptographicMessageSyntax](#) ()
- bool [Decrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*decrypt content from a CMS envelopedData structure*
- bool [Encrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*create a CMS envelopedData structure*
- [CipherTypes](#) [GetCipherType](#) () const
- bool [GetInitialized](#) () const
- bool [ParseCertificateFile](#) (const char \*filename)
- bool [ParseKeyFile](#) (const char \*filename)
- void [SetCipherType](#) ([CipherTypes](#) type)
- bool [SetPassword](#) (const char \*pass, size\_t passLen)

## Additional Inherited Members

### 27.41.1 Constructor & Destructor Documentation

27.41.1.1 `gdcM::CAPICryptographicMessageSyntax::CAPICryptographicMessageSyntax ( )`

27.41.1.2 `gdcM::CAPICryptographicMessageSyntax::~~CAPICryptographicMessageSyntax ( )`

### 27.41.2 Member Function Documentation

27.41.2.1 `bool gdcM::CAPICryptographicMessageSyntax::Decrypt ( char * output, size_t & outlen, const char * array, size_t len ) const [virtual]`

decrypt content from a CMS envelopedData structure

Implements [gdcM::CryptographicMessageSyntax](#).

27.41.2.2 `bool gdcM::CAPICryptographicMessageSyntax::Encrypt ( char * output, size_t & outlen, const char * array, size_t len ) const [virtual]`

create a CMS envelopedData structure

Implements [gdcM::CryptographicMessageSyntax](#).

27.41.2.3 `CipherTypes gdcM::CAPICryptographicMessageSyntax::GetCipherType ( ) const [virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.41.2.4 `bool gdcM::CAPICryptographicMessageSyntax::GetInitialized ( ) const [inline]`

27.41.2.5 `bool gdcM::CAPICryptographicMessageSyntax::ParseCertificateFile ( const char * filename ) [virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.41.2.6 `bool gdcM::CAPICryptographicMessageSyntax::ParseKeyFile ( const char * filename ) [virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.41.2.7 `void gdcM::CAPICryptographicMessageSyntax::SetCipherType ( CipherTypes type )`

27.41.2.8 `bool gdcM::CAPICryptographicMessageSyntax::SetPassword ( const char * pass, size_t passLen ) [virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

The documentation for this class was generated from the following file:


- [gdcMCAPICryptographicMessageSyntax.h](#)

## 27.42 gdcm::network::CEchoRQ Class Reference


[CEchoRQ](#) this file defines the messages for the cecho action.

```
#include <gdcmCEchoMessages.h>
```

Inheritance diagram for gdcm::network::CEchoRQ:

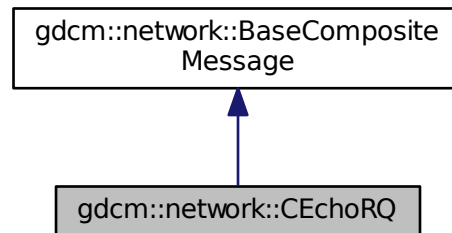

Collaboration diagram for gdcm::network::CEchoRQ:

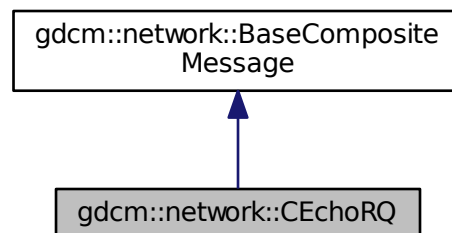

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDV (const ULConnection &inConnection, const BaseRootQuery *inRootQuery)`

### Public Attributes

- `UIComp AffectedSOPClassUID`
- `uint16_t MessageID`

### 27.42.1 Detailed Description

[CEchoRQ](#) this file defines the messages for the cecho action.

### 27.42.2 Member Function Documentation

27.42.2.1 `std::vector<PresentationDataValue> gdcmm::network::CEchoRQ::ConstructPDV ( const ULConnection & inConnection, const BaseRootQuery * inRootQuery ) [virtual]`

Implements [gdcmm::network::BaseCompositeMessage](#).

### 27.42.3 Member Data Documentation

27.42.3.1 `UIComp gdcmm::network::CEchoRQ::AffectedSOPClassUID`

27.42.3.2 `uint16_t gdcmm::network::CEchoRQ::MessageID`

The documentation for this class was generated from the following files:


- [gdcmmCEchoMessages.h](#)
- [gdcmmDIMSE.h](#)

## 27.43 gdcmm::network::CEchoRSP Class Reference


[CEchoRSP](#) this file defines the messages for the cecho action.

```
#include <gdcmmCEchoMessages.h>
```

Inheritance diagram for `gdcmm::network::CEchoRSP`:

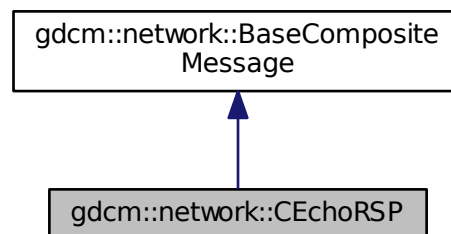

Collaboration diagram for gdcmm::network::CEchoRSP:

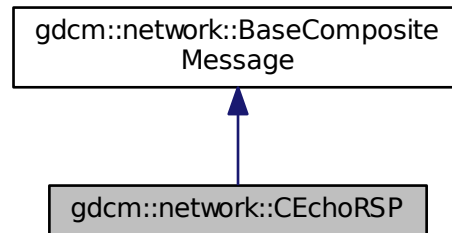

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDVByDataSet (const DataSet *inDataSet)`

#### 27.43.1 Detailed Description

[CEchoRSP](#) this file defines the messages for the cecho action.

#### 27.43.2 Member Function Documentation

27.43.2.1 `std::vector<PresentationDataValue> gdcmm::network::CEchoRSP::ConstructPDVByDataSet ( const DataSet *inDataSet )`

The documentation for this class was generated from the following file:

- [gdcmmCEchoMessages.h](#)

## 27.44 gdcmm::network::CFind Class Reference

```
#include <gdcmmDIMSE.h>
```

#### 27.44.1 Detailed Description

PS 3.4 - 2009 [Table B.2-1](#) C-STORE STATUS

The documentation for this class was generated from the following file:


- [gdcmmDIMSE.h](#)

## 27.45 gdcmm::network::CFindCancelRQ Class Reference


[CFindCancelRQ](#) this file defines the messages for the cfind action.

```
#include <gdcmmCFindMessages.h>
```

Inheritance diagram for gdcmm::network::CFindCancelRQ:

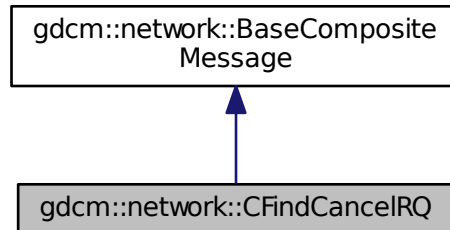

Collaboration diagram for gdcmm::network::CFindCancelRQ:

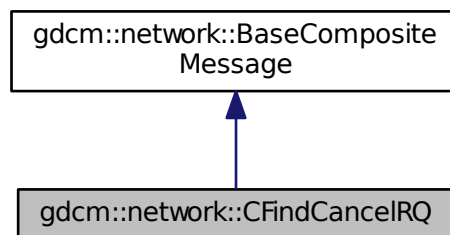

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDVByDataSet` (const [DataSet](#) \*inDataSet)

### 27.45.1 Detailed Description

[CFindCancelRQ](#) this file defines the messages for the cfind action.

### 27.45.2 Member Function Documentation



27.45.2.1 `std::vector<PresentationDataValue> gdcm::network::CFindCancelRQ::ConstructPDVByDataSet ( const DataSet * inDataSet )`

The documentation for this class was generated from the following file:


- [gdcmCFindMessages.h](#)

## 27.46 gdcm::network::CFindRQ Class Reference


[CFindRQ](#) this file defines the messages for the cfind action.

```
#include <gdcmCFindMessages.h>
```

Inheritance diagram for `gdcm::network::CFindRQ`:

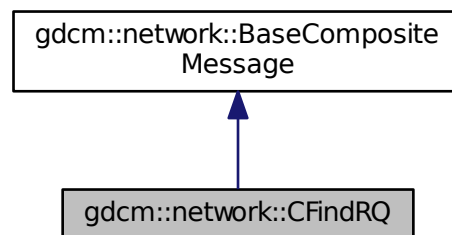

Collaboration diagram for `gdcm::network::CFindRQ`:

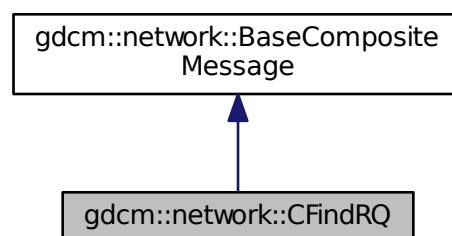

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDV (const ULConnection &inConnection, const BaseRootQuery *inRootQuery)`

### 27.46.1 Detailed Description

[CFindRQ](#) this file defines the messages for the cfind action.

### 27.46.2 Member Function Documentation

27.46.2.1 `std::vector<PresentationDataValue> gdcmm::network::CFindRQ::ConstructPDV ( const ULConnection & inConnection, const BaseRootQuery * inRootQuery ) [virtual]`

Implements [gdcmm::network::BaseCompositeMessage](#).

The documentation for this class was generated from the following file:


- [gdcmmCFindMessages.h](#)

## 27.47 gdcmm::network::CFindRSP Class Reference


[CFindRSP](#) this file defines the messages for the cfind action.

```
#include <gdcmmCFindMessages.h>
```

Inheritance diagram for gdcmm::network::CFindRSP:

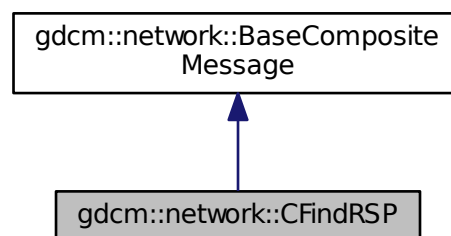

Collaboration diagram for gdcm::network::CFindRSP:

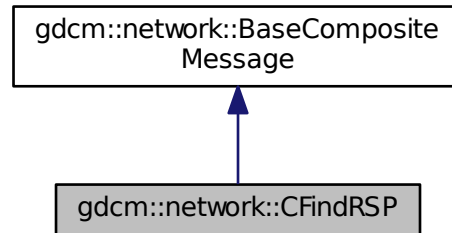

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDVByDataSet (const DataSet *inDataSet)`

#### 27.47.1 Detailed Description

[CFindRSP](#) this file defines the messages for the cfind action.

#### 27.47.2 Member Function Documentation

27.47.2.1 `std::vector<PresentationDataValue> gdcm::network::CFindRSP::ConstructPDVByDataSet ( const DataSet *inDataSet )`


The documentation for this class was generated from the following file:

- [gdcmCFindMessages.h](#)


## 27.48 gdcm::network::CMoveCancelRq Class Reference

```
#include <gdcmCMoveMessages.h>
```

Inheritance diagram for `gdcm::network::CMoveCancelRq`:

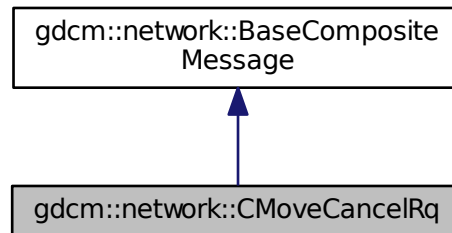

Collaboration diagram for `gdcm::network::CMoveCancelRq`:

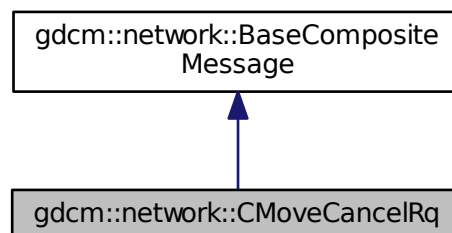

## Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDVByDataSet (const DataSet *inDataSet)`

### 27.48.1 Member Function Documentation

27.48.1.1 `std::vector<PresentationDataValue> gdcm::network::CMoveCancelRq::ConstructPDVByDataSet ( const DataSet *inDataSet )`

The documentation for this class was generated from the following file:


- [gdcmCMoveMessages.h](#)

## 27.49 gdcm::network::CMoveRQ Class Reference


[CMoveRQ](#) this file defines the messages for the cmove action.

```
#include <gdcmCMoveMessages.h>
```

Inheritance diagram for gdcm::network::CMoveRQ:

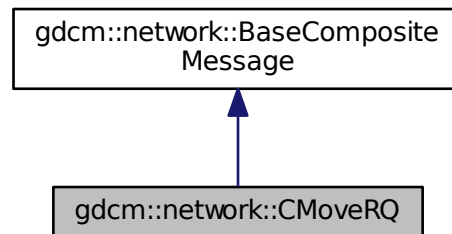

Collaboration diagram for gdcm::network::CMoveRQ:

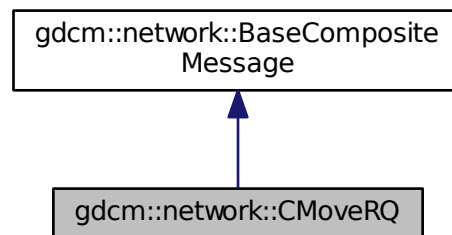

### Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDV` (const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \*inRootQuery)

### 27.49.1 Detailed Description

[CMoveRQ](#) this file defines the messages for the cmove action.

### 27.49.2 Member Function Documentation

27.49.2.1 `std::vector<PresentationDataValue> gdcm::network::CMoveRQ::ConstructPDV ( const ULConnection & inConnection, const BaseRootQuery * inRootQuery ) [virtual]`

Implements [gdcm::network::BaseCompositeMessage](#).

The documentation for this class was generated from the following file:


- [gdcmCMoveMessages.h](#)

## 27.50 gdcm::network::CMoveRSP Class Reference


[CMoveRSP](#) this file defines the messages for the cmove action.

```
#include <gdcmCMoveMessages.h>
```

Inheritance diagram for `gdcm::network::CMoveRSP`:

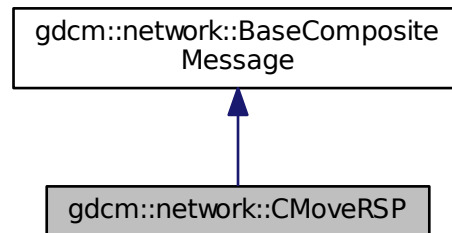

Collaboration diagram for `gdcm::network::CMoveRSP`:

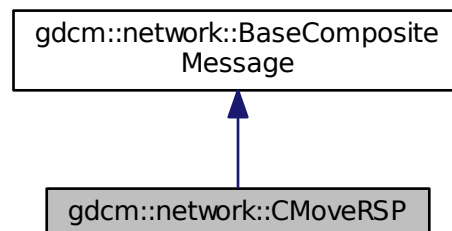

## Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDVByDataSet (const DataSet *inDataSet)`

### 27.50.1 Detailed Description

[CMoveRSP](#) this file defines the messages for the cmove action.

### 27.50.2 Member Function Documentation

27.50.2.1 `std::vector<PresentationDataValue> gdcm::network::CMoveRSP::ConstructPDVByDataSet ( const DataSet *inDataSet )`

The documentation for this class was generated from the following file:


- [gdcmCMoveMessages.h](#)

## 27.51 gdcm::Codec Class Reference


[Codec](#) class.

```
#include <gdcmCodec.h>
```

Inheritance diagram for `gdcm::Codec`:

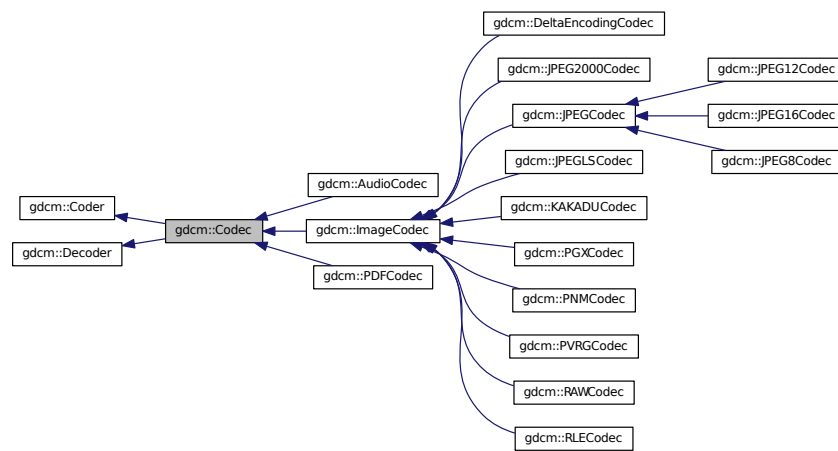

Collaboration diagram for `gdcm::Codec`:

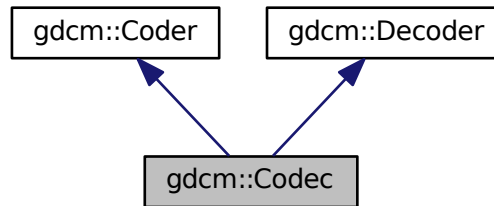

## Additional Inherited Members

### 27.51.1 Detailed Description

`Codec` class.

The documentation for this class was generated from the following file:


- [gdcmCodec.h](#)

## 27.52 `gdcm::Coder` Class Reference

`Coder`.

```
#include <gdcmCoder.h>
```

Inheritance diagram for `gdcm::Coder`:

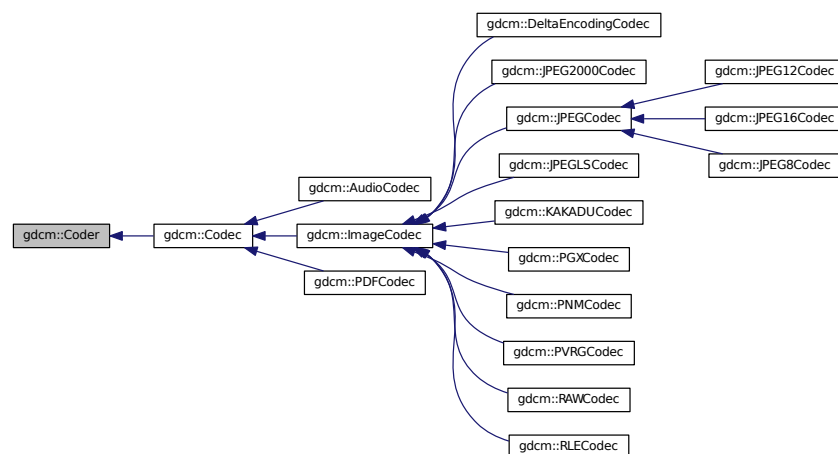



## Public Member Functions

- virtual [~Coder](#) ()
- virtual bool [CanCode](#) ([TransferSyntax](#) const &) const =0  
*Return whether this coder support this transfer syntax (can code it)*
- virtual bool [Code](#) ([DataElement](#) const &in\_, [DataElement](#) &out\_)  
*Code.*

## Protected Member Functions

- virtual bool [InternalCode](#) (const char \*bv, unsigned long len, std::ostream &os)

### 27.52.1 Detailed Description

[Coder](#).

### 27.52.2 Constructor & Destructor Documentation

27.52.2.1 virtual gdcm::Coder::~Coder ( ) [inline], [virtual]

### 27.52.3 Member Function Documentation

27.52.3.1 virtual bool gdcm::Coder::CanCode ( [TransferSyntax](#) const & ) const [pure virtual]

Return whether this coder support this transfer syntax (can code it)

Implemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::JPEGLSCodec](#), [gdcm::ImageCodec](#), [gdcm::PNMCodec](#), [gdcm::PGXCodec](#), [gdcm::KAKADUCodec](#), [gdcm::RAWCodec](#), [gdcm::AudioCodec](#), and [gdcm::PDFCodec](#).

27.52.3.2 virtual bool gdcm::Coder::Code ( [DataElement](#) const & in\_, [DataElement](#) & out\_ ) [inline], [virtual]

Code.

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::JPEGLSCodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::KAKADUCodec](#), and [gdcm::RAWCodec](#).

27.52.3.3 virtual bool gdcm::Coder::InternalCode ( const char \* bv, unsigned long len, std::ostream & os ) [inline], [protected], [virtual]

Reimplemented in [gdcm::JPEG12Codec](#), [gdcm::JPEG16Codec](#), and [gdcm::JPEG8Codec](#).

The documentation for this class was generated from the following file:

- [gdcmCoder.h](#)

## 27.53 gdcm::CodeString Class Reference

[CodeString](#) This is an implementation of DICOM [VR](#): CS The ctor will properly Trim so that operator== is correct.

```
#include <gdcmCodeString.h>
```

## Public Types

- typedef [InternalClass::const\\_iterator](#) const\_iterator
- typedef [InternalClass::const\\_reference](#) const\_reference
- typedef [InternalClass::const\\_reverse\\_iterator](#) const\_reverse\_iterator
- typedef [InternalClass::difference\\_type](#) difference\_type
- typedef [InternalClass::iterator](#) iterator
- typedef [InternalClass::pointer](#) pointer
- typedef [InternalClass::reference](#) reference
- typedef [InternalClass::reverse\\_iterator](#) reverse\_iterator
- typedef [InternalClass::size\\_type](#) size\_type
- typedef [InternalClass::value\\_type](#) value\_type

## Public Member Functions

- [CodeString](#) ()  
*[CodeString](#) constructors.*
- [CodeString](#) (const [value\\_type](#) \*s)
- [CodeString](#) (const [value\\_type](#) \*s, [size\\_type](#) n)
- [CodeString](#) (const [InternalClass](#) &s, [size\\_type](#) pos=0, [size\\_type](#) n=[InternalClass::npos](#))
- std::string [GetAsString](#) () const  
*Return the full code string as std::string.*
- bool [IsValid](#) () const  
*Check if [CodeString](#) obj is correct..*
- [size\\_type](#) [Size](#) () const  
*Return the size of the string.*

## Protected Member Functions

- std::string [TrimInternal](#) () const

## Friends

- bool [operator!=](#) (const [CodeString](#) &ref, const [CodeString](#) &cs)
- std::ostream & [operator<<](#) (std::ostream &os, const [CodeString](#) &str)
- bool [operator==](#) (const [CodeString](#) &ref, const [CodeString](#) &cs)

### 27.53.1 Detailed Description

[CodeString](#) This is an implementation of DICOM [VR: CS](#) The ctor will properly Trim so that operator== is correct.

**Note**

the ctor of [CodeString](#) will Trim the string on the fly so as to remove the extra leading and ending spaces. However it will not perform validation on the fly ([CodeString](#) obj can contains invalid char such as lower cases). This design was chosen to be a little tolerant to broken DICOM implementation, and thus allow user to compare lower case CS from there input file without the need to first rewrite them to get rid of invalid character (validation is a different operation from searching, querying).

**Warning**

when writing out DICOM file it is highly recommended to perform the [IsValid\(\)](#) call, at least to check that the length of the string match the definition in the standard.

**27.53.2 Member Typedef Documentation**

27.53.2.1 `typedef InternalClass::const_iterator gdcm::CodeString::const_iterator`

27.53.2.2 `typedef InternalClass::const_reference gdcm::CodeString::const_reference`

27.53.2.3 `typedef InternalClass::const_reverse_iterator gdcm::CodeString::const_reverse_iterator`

27.53.2.4 `typedef InternalClass::difference_type gdcm::CodeString::difference_type`

27.53.2.5 `typedef InternalClass::iterator gdcm::CodeString::iterator`

27.53.2.6 `typedef InternalClass::pointer gdcm::CodeString::pointer`

27.53.2.7 `typedef InternalClass::reference gdcm::CodeString::reference`

27.53.2.8 `typedef InternalClass::reverse_iterator gdcm::CodeString::reverse_iterator`

27.53.2.9 `typedef InternalClass::size_type gdcm::CodeString::size_type`

27.53.2.10 `typedef InternalClass::value_type gdcm::CodeString::value_type`

**27.53.3 Constructor & Destructor Documentation**

27.53.3.1 `gdcm::CodeString::CodeString ( ) [inline]`

[CodeString](#) constructors.

27.53.3.2 `gdcm::CodeString::CodeString ( const value_type * s ) [inline]`

27.53.3.3 `gdcm::CodeString::CodeString ( const value_type * s, size_type n ) [inline]`

27.53.3.4 `gdcm::CodeString::CodeString ( const InternalClass & s, size_type pos = 0, size_type n = InternalClass::npos ) [inline]`

**27.53.4 Member Function Documentation**

27.53.4.1 `std::string gdcM::CodeString::GetAsString ( ) const` `[inline]`

Return the full code string as `std::string`.

27.53.4.2 `bool gdcM::CodeString::IsValid ( ) const`

Check if [CodeString](#) obj is correct..

27.53.4.3 `size_type gdcM::CodeString::Size ( ) const` `[inline]`

Return the size of the string.

27.53.4.4 `std::string gdcM::CodeString::TrimInternal ( ) const` `[inline],[protected]`

## 27.53.5 Friends And Related Function Documentation

27.53.5.1 `bool operator!= ( const CodeString & ref, const CodeString & cs )` `[friend]`

27.53.5.2 `std::ostream& operator<< ( std::ostream & os, const CodeString & str )` `[friend]`

27.53.5.3 `bool operator== ( const CodeString & ref, const CodeString & cs )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcMCodeString.h](#)

## 27.54 gdcM::Command Class Reference


[Command](#) superclass for callback/observer methods.

```
#include <gdcMCommand.h>
```

Inheritance diagram for gdcmm::Command:

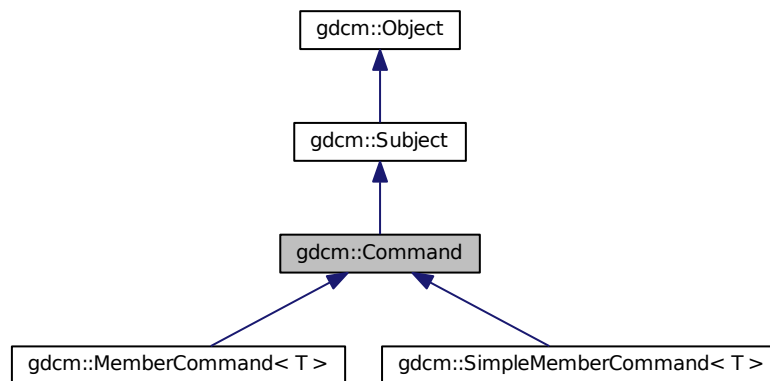

Collaboration diagram for gdcmm::Command:

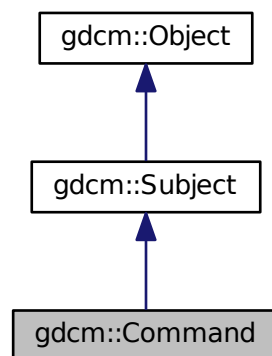

## Public Member Functions

- virtual void `Execute` (`Subject` \*caller, const `Event` &event)=0  
*Abstract method that defines the action to be taken by the command.*
- virtual void `Execute` (const `Subject` \*caller, const `Event` &event)=0

## Protected Member Functions

- `Command` ()
- `~Command` ()

### 27.54.1 Detailed Description

[Command](#) superclass for callback/observer methods.

See also

[Subject](#)

### 27.54.2 Constructor & Destructor Documentation

27.54.2.1 `gdcm::Command::Command ( )` [protected]

27.54.2.2 `gdcm::Command::~~Command ( )` [protected]

### 27.54.3 Member Function Documentation

27.54.3.1 `virtual void gdcm::Command::Execute ( Subject * caller, const Event & event )` [pure virtual]

Abstract method that defines the action to be taken by the command.

Implemented in [gdcm::SimpleMemberCommand< T >](#), and [gdcm::MemberCommand< T >](#).

27.54.3.2 `virtual void gdcm::Command::Execute ( const Subject * caller, const Event & event )` [pure virtual]

Abstract method that defines the action to be taken by the command. This variant is expected to be used when requests comes from a const [Object](#)

Implemented in [gdcm::SimpleMemberCommand< T >](#), and [gdcm::MemberCommand< T >](#).

The documentation for this class was generated from the following file:


- [gdcmCommand.h](#)

## 27.55 gdcm::CommandDataSet Class Reference


Class to represent a [Command DataSet](#).

```
#include <gdcmCommandDataSet.h>
```

Inheritance diagram for gdcm::CommandDataSet:

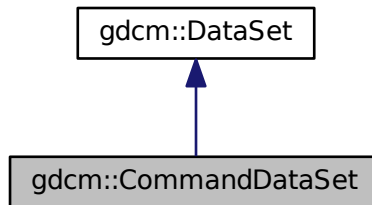

Collaboration diagram for gdcm::CommandDataSet:

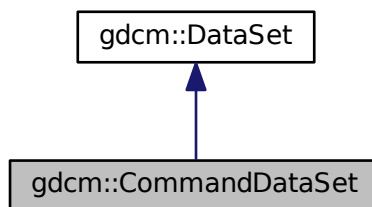

## Public Member Functions

- [CommandDataSet](#) ()
- [~CommandDataSet](#) ()
- void [Insert](#) (const [DataElement](#) &de)
- std::istream & [Read](#) (std::istream &is)  
*Read.*
- void [Replace](#) (const [DataElement](#) &de)
- std::ostream & [Write](#) (std::ostream &os) const  
*Write.*

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [CommandDataSet](#) &\_val)

## Additional Inherited Members

### 27.55.1 Detailed Description

Class to represent a [Command DataSet](#).

See also

[DataSet](#)

### 27.55.2 Constructor & Destructor Documentation

27.55.2.1 `gdcm::CommandDataSet::CommandDataSet ( )` `[inline]`

27.55.2.2 `gdcm::CommandDataSet::~~CommandDataSet ( )` `[inline]`

### 27.55.3 Member Function Documentation

27.55.3.1 `void gdcm::CommandDataSet::Insert ( const DataElement & de )` `[inline]`

References `gdcmErrorMacro`, `gdcm::Tag::GetGroup()`, and `gdcm::DataElement::GetTag()`.

27.55.3.2 `std::istream& gdcm::CommandDataSet::Read ( std::istream & is )`

Read.

27.55.3.3 `void gdcm::CommandDataSet::Replace ( const DataElement & de )` `[inline]`

References `gdcm::DataElement::GetTag()`.

27.55.3.4 `std::ostream& gdcm::CommandDataSet::Write ( std::ostream & os ) const`

Write.

### 27.55.4 Friends And Related Function Documentation

27.55.4.1 `std::ostream& operator<< ( std::ostream & _os, const CommandDataSet & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmCommandDataSet.h](#)

## 27.56 gdcm::network::CompositeMessageFactory Class Reference

[CompositeMessageFactory](#) This class constructs PDataPDUs, but that have been specifically constructed for the composite DICOM services (C-Echo, C-Find, C-Get, C-Move, and C-Store). It will also handle parsing the incoming data to determine which of the CompositePDUs the incoming data is, and so therefore allowing the scu to determine what to do with incoming data (if acting as a storescp server, for instance).



```
#include <gdcmCompositeMessageFactory.h>
```

### Static Public Member Functions

- static std::vector< [PresentationDataValue](#) > [ConstructCEchoRQ](#) (const [ULConnection](#) &inConnection)
- static std::vector< [PresentationDataValue](#) > [ConstructCFindRQ](#) (const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \*inRootQuery)
- static std::vector< [PresentationDataValue](#) > [ConstructCMoveRQ](#) (const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \*inRootQuery)
- static std::vector< [PresentationDataValue](#) > [ConstructCStoreRQ](#) (const [ULConnection](#) &inConnection, const [File](#) &file)
- static std::vector< [PresentationDataValue](#) > [ConstructCStoreRSP](#) (const [DataSet](#) \*inDataSet, const [BasePDU](#) \*inPC)

#### 27.56.1 Detailed Description

[CompositeMessageFactory](#) This class constructs PDataPDUs, but that have been specifically constructed for the composite DICOM services (C-Echo, C-Find, C-Get, C-Move, and C-Store). It will also handle parsing the incoming data to determine which of the CompositePDUs the incoming data is, and so therefore allowing the scu to determine what to do with incoming data (if acting as a storescp server, for instance).

#### 27.56.2 Member Function Documentation

27.56.2.1 static std::vector<[PresentationDataValue](#)> gdcm::network::CompositeMessageFactory::ConstructCEchoRQ ( const [ULConnection](#) & *inConnection* ) [static]

27.56.2.2 static std::vector<[PresentationDataValue](#)> gdcm::network::CompositeMessageFactory::ConstructCFindRQ ( const [ULConnection](#) & *inConnection*, const [BaseRootQuery](#) \* *inRootQuery* ) [static]

27.56.2.3 static std::vector<[PresentationDataValue](#)> gdcm::network::CompositeMessageFactory::ConstructCMoveRQ ( const [ULConnection](#) & *inConnection*, const [BaseRootQuery](#) \* *inRootQuery* ) [static]

27.56.2.4 static std::vector<[PresentationDataValue](#)> gdcm::network::CompositeMessageFactory::ConstructCStoreRQ ( const [ULConnection](#) & *inConnection*, const [File](#) & *file* ) [static]

27.56.2.5 static std::vector<[PresentationDataValue](#)> gdcm::network::CompositeMessageFactory::ConstructCStoreRSP ( const [DataSet](#) \* *inDataSet*, const [BasePDU](#) \* *inPC* ) [static]

The documentation for this class was generated from the following file:

- [gdcmCompositeMessageFactory.h](#)

## 27.57 gdcm::CompositeNetworkFunctions Class Reference

**Composite Network Functions** These functions provide a generic API to the DICOM functions implemented in GDCM. Advanced users can use this code as a template for building their own versions of these functions (for instance, to provide progress bars or some other way of handling returned query information), but for most users, these functions should be sufficient to interface with a PACS to a local machine. Note that these functions are not contained within a static class or some other class-style interface, because multiple connections can be instantiated in the same program.

The DICOM standard is much more function oriented rather than class oriented in this instance, so the design of this API reflects that functional approach. These functions implements the following SCU operations:

```
#include <gdcmCompositeNetworkFunctions.h>
```

## Public Types

- typedef std::vector< [KeyValuePairType](#) > [KeyValuePairArrayType](#)
- typedef std::pair< [Tag](#), std::string > [KeyValuePairType](#)

## Static Public Member Functions

- static bool [CEcho](#) (const char \*remote, uint16\_t portno, const char \*aetitle=NULL, const char \*call=NULL)
- static bool [CFind](#) (const char \*remote, uint16\_t portno, const [BaseRootQuery](#) \*query, std::vector< [DataSet](#) > &retDataSets, const char \*aetitle=NULL, const char \*call=NULL)
- static bool [CMove](#) (const char \*remote, uint16\_t portno, const [BaseRootQuery](#) \*query, uint16\_t portscp, const char \*aetitle=NULL, const char \*call=NULL, const char \*outputdir=NULL)
- static [BaseRootQuery](#) \* [ConstructQuery](#) ([ERootType](#) inRootType, [EQueryLevel](#) inQueryLevel, const [DataSet](#) &queryds, bool inMove=false)
- static [BaseRootQuery](#) \* [ConstructQuery](#) ([ERootType](#) inRootType, [EQueryLevel](#) inQueryLevel, const [KeyValuePairArrayType](#) &keys, bool inMove=false)
- static bool [CStore](#) (const char \*remote, uint16\_t portno, const [Directory::FileNamesType](#) &filenames, const char \*aetitle=NULL, const char \*call=NULL)

### 27.57.1 Detailed Description

Composite Network Functions These functions provide a generic API to the DICOM functions implemented in GDCM. Advanced users can use this code as a template for building their own versions of these functions (for instance, to provide progress bars or some other way of handling returned query information), but for most users, these functions should be sufficient to interface with a PACS to a local machine. Note that these functions are not contained within a static class or some other class-style interface, because multiple connections can be instantiated in the same program. The DICOM standard is much more function oriented rather than class oriented in this instance, so the design of this API reflects that functional approach. These functions implements the following SCU operations:

- C-ECHO SCU
- C-FIND SCU
- C-STORE SCU
- C-MOVE SCU (+internal C-STORE SCP)

### 27.57.2 Member Typedef Documentation

27.57.2.1 typedef std::vector< [KeyValuePairType](#) > [gdcm::CompositeNetworkFunctions::KeyValuePairArrayType](#)

27.57.2.2 typedef std::pair< [Tag](#), std::string > [gdcm::CompositeNetworkFunctions::KeyValuePairType](#)

### 27.57.3 Member Function Documentation

27.57.3.1 `static bool gdcm::CompositeNetworkFunctions::CEcho ( const char * remote, uint16_t portno, const char * aetitle = NULL, const char * call = NULL ) [static]`

The most basic network function. Use this function to ensure that the remote server is responding on the given IP and port number as expected.

## Parameters

<i>aetitle</i>	when not set will default to 'GDCMSCU'
<i>call</i>	when not set will default to 'ANY-SCP'

## Warning

This is an error to set remote to NULL or portno to 0

## Returns

true if it worked.

**27.57.3.2** `static bool gdcmm::CompositeNetworkFunctions::CFind ( const char * remote, uint16_t portno, const BaseRootQuery * query, std::vector< DataSet > & retDataSets, const char * aetitle = NULL, const char * call = NULL ) [static]`

This function will use the provided query to determine what files a remote server contains that match the query strings. The return is a vector of datasets that contain tags as reported by the server. If the dataset is empty, then it is possible that an error condition was encountered; in which case, the user should monitor the error and warning streams.

## Parameters

<i>aetitle</i>	when not set will default to 'GDCMSCU'
<i>call</i>	when not set will default to 'ANY-SCP'

## Warning

This is an error to set remote to NULL or portno to 0

## Returns

true if it worked.

**27.57.3.3** `static bool gdcmm::CompositeNetworkFunctions::CMove ( const char * remote, uint16_t portno, const BaseRootQuery * query, uint16_t portscp, const char * aetitle = NULL, const char * call = NULL, const char * outputdir = NULL ) [static]`

This function will use the provided query to get files from a remote server. NOTE that this functionality is essentially equivalent to C-GET in the DICOM standard; however, C-GET has been deprecated, so this function allows for the user to ask a remote server for files matching a query and return them to the local machine. Files will be written to the given output directory. If the operation succeeds, the function returns true. This function is a prime candidate for being overwritten by expert users; if the datasets should remain in memory, for instance, that behavior can be changed by creating a user-level version of this function.

## Parameters

<i>aetitle</i>	when not set will default to 'GDCMSCU'
<i>call</i>	when not set will default to 'ANY-SCP' This is an error to set remote to NULL or portno to 0 when

<i>outputdir</i>	is not set default to current dir ('.')
------------------	---

**Returns**

true if it worked.

**27.57.3.4** static **BaseRootQuery\*** gdcm::CompositeNetworkFunctions::ConstructQuery ( **ERootType** *inRootType*, **EQueryLevel** *inQueryLevel*, const **DataSet** & *queryds*, bool *inMove* = false ) [static]

This function will take a list of strings and tags and fill in a query that can be used for either CFind or CMove (depending on the input boolean

**Parameters**

<i>inMove</i> ).	Note that the caller is responsible for deleting the constructed query. This function is used to build both a move and a find query (true for inMove if it's move, false if it's find)
------------------	--

**27.57.3.5** static **BaseRootQuery\*** gdcm::CompositeNetworkFunctions::ConstructQuery ( **ERootType** *inRootType*, **EQueryLevel** *inQueryLevel*, const **KeyValuePairArrayType** & *keys*, bool *inMove* = false ) [static]

**Deprecated**

**27.57.3.6** static bool gdcm::CompositeNetworkFunctions::CStore ( const char \* *remote*, uint16\_t *portno*, const **Directory::FileNamesType** & *filenames*, const char \* *aetitle* = NULL, const char \* *call* = NULL ) [static]

This function will place the provided files into the remote server. The function returns true if it worked for all files.

**Warning**

the server side can refuse an association on a given file

**Parameters**

<i>aetitle</i>	when not set will default to 'GDCMSCU'
<i>call</i>	when not set will default to 'ANY-SCP'

**Warning**

This is an error to set remote to NULL or portno to 0

**Returns**

true if it worked for all files

The documentation for this class was generated from the following file:

- [gdcmCompositeNetworkFunctions.h](#)

## 27.58 gdcm::ConstCharWrapper Class Reference

Do not use me.

```
#include <gdcmConstCharWrapper.h>
```

### Public Member Functions

- [ConstCharWrapper](#) (const char \*i=0)
- [operator const char \\* \(\)](#) const

### 27.58.1 Detailed Description

Do not use me.

### 27.58.2 Constructor & Destructor Documentation

27.58.2.1 `gdcm::ConstCharWrapper::ConstCharWrapper ( const char * i = 0 )` `[inline]`

### 27.58.3 Member Function Documentation

27.58.3.1 `gdcm::ConstCharWrapper::operator const char * ( )` const `[inline]`

The documentation for this class was generated from the following file:


- [gdcmConstCharWrapper.h](#)

## 27.59 gdcm::CP246ExplicitDataElement Class Reference


Class to read/write a [DataElement](#) as CP246Explicit Data [Element](#).

```
#include <gdcmCP246ExplicitDataElement.h>
```

Inheritance diagram for `gdcm::CP246ExplicitDataElement`:

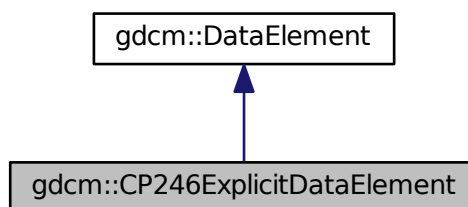

Collaboration diagram for gdcm::CP246ExplicitDataElement:

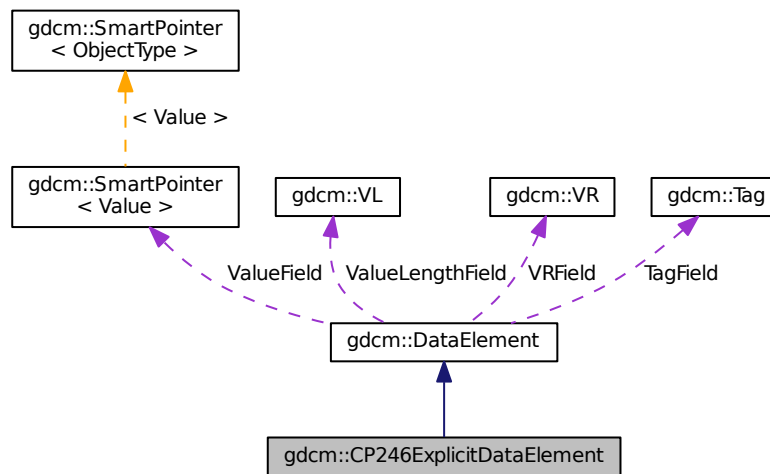

## Public Member Functions

- [VL GetLength](#) () const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadPreValue](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadValue](#) (std::istream &is, bool readvalues=true)
- template<typename TSwap >  
std::istream & [ReadWithLength](#) (std::istream &is, [VL](#) &length)

## Additional Inherited Members

### 27.59.1 Detailed Description

Class to read/write a [DataElement](#) as CP246Explicit Data [Element](#).

#### Note

Some system are producing SQ, declare them as UN, but encode the SQ as 'Explicit' instead of Implicit

### 27.59.2 Member Function Documentation

27.59.2.1 [VL](#) gdcm::CP246ExplicitDataElement::GetLength ( ) const

27.59.2.2 template<typename TSwap > std::istream& gdcm::CP246ExplicitDataElement::Read ( std::istream & is )

27.59.2.3 `template<typename TSwap> std::istream& gdcM::CP246ExplicitDataElement::ReadPreValue ( std::istream & is )`

27.59.2.4 `template<typename TSwap> std::istream& gdcM::CP246ExplicitDataElement::ReadValue ( std::istream & is, bool readvalues = true )`

27.59.2.5 `template<typename TSwap> std::istream& gdcM::CP246ExplicitDataElement::ReadWithLength ( std::istream & is, VL & length )`

The documentation for this class was generated from the following file:


- [gdcMCP246ExplicitDataElement.h](#)

## 27.60 gdcM::CryptoFactory Class Reference

Class to do handle the crypto factory.

```
#include <gdcMCryptoFactory.h>
```

Inheritance diagram for gdcM::CryptoFactory:

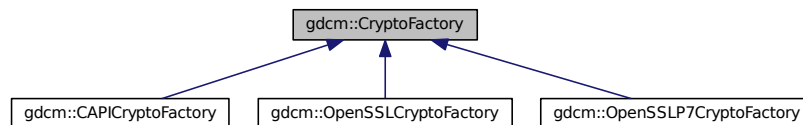

### Public Types

- enum `CryptoLib` {  
`DEFAULT` = 0,  
`OPENSSL` = 1,  
`CAPI` = 2,  
`OPENSSL7` = 3 }

### Public Member Functions

- virtual `CryptographicMessageSyntax * CreateCMSProvider ()`=0

### Static Public Member Functions

- static `CryptoFactory * GetFactoryInstance (CryptoLib id=DEFAULT)`

### Protected Member Functions

- `CryptoFactory (CryptoLib id)`
- `CryptoFactory ()`
- `~CryptoFactory ()`



### 27.60.1 Detailed Description

Class to do handle the crypto factory.

GDCM needs to access in a platform independant way the user specified crypto engine. It can be:

- CAPI (windows only)
- OPENSSL (portable)
- OPENSSL7 (portable) By default the factory will try: CAPI if on windows OPENSSL if possible OPENSSL7 when older OpenSSL is used.

### 27.60.2 Member Enumeration Documentation

#### 27.60.2.1 enum gdcm::CryptoFactory::CryptoLib

Enumerator

***DEFAULT***

***OPENSSL***

***CAPI***

***OPENSSL7***

### 27.60.3 Constructor & Destructor Documentation

27.60.3.1 `gdcm::CryptoFactory::CryptoFactory ( CryptoLib id )` `[inline]`, `[protected]`

27.60.3.2 `gdcm::CryptoFactory::CryptoFactory ( )` `[inline]`, `[protected]`

27.60.3.3 `gdcm::CryptoFactory::~~CryptoFactory ( )` `[inline]`, `[protected]`

### 27.60.4 Member Function Documentation

27.60.4.1 `virtual CryptographicMessageSyntax* gdcm::CryptoFactory::CreateCMSProvider ( )` `[pure virtual]`

Implemented in [gdcm::OpenSSLCryptoFactory](#), [gdcm::OpenSSL7CryptoFactory](#), and [gdcm::CAPICryptoFactory](#).

27.60.4.2 `static CryptoFactory* gdcm::CryptoFactory::GetFactoryInstance ( CryptoLib id = DEFAULT )` `[static]`

Examples:

[BasicAnonymizer.cs](#), and [ClinicalTrialIdentificationWorkflow.cs](#).


The documentation for this class was generated from the following file:

- [gdcmCryptoFactory.h](#)

## 27.61 gdcmm::CryptographicMessageSyntax Class Reference

```
#include <gdcmmCryptographicMessageSyntax.h>
```

Inheritance diagram for gdcmm::CryptographicMessageSyntax:

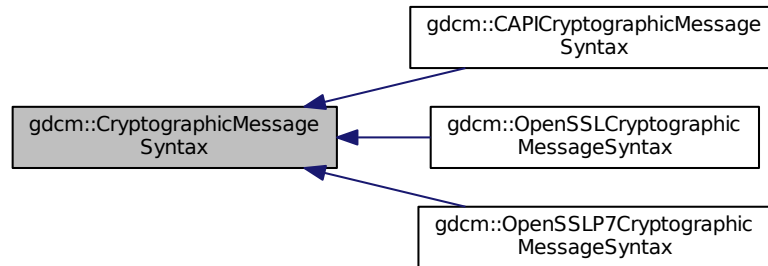

### Public Types

- enum [CipherTypes](#) {  
[DES3\\_CIPHER](#),  
[AES128\\_CIPHER](#),  
[AES192\\_CIPHER](#),  
[AES256\\_CIPHER](#) }

### Public Member Functions

- [CryptographicMessageSyntax](#) ()
- virtual [~CryptographicMessageSyntax](#) ()
- virtual bool [Decrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const =0  
*decrypt content from a CMS envelopedData structure*
- virtual bool [Encrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const =0  
*create a CMS envelopedData structure*
- virtual [CipherTypes GetCipherType](#) () const =0
- virtual bool [ParseCertificateFile](#) (const char \*filename)=0
- virtual bool [ParseKeyFile](#) (const char \*filename)=0
- virtual void [SetCipherType](#) ([CipherTypes](#) type)=0
- virtual bool [SetPassword](#) (const char \*pass, size\_t passLen)=0

### 27.61.1 Member Enumeration Documentation

#### 27.61.1.1 enum gdcmm::CryptographicMessageSyntax::CipherTypes

Enumerator

***DES3\_CIPHER***  
***AES128\_CIPHER***

***AES192\_CIPHER******AES256\_CIPHER***

## 27.61.2 Constructor & Destructor Documentation

27.61.2.1 `gdcmm::CryptographicMessageSyntax::CryptographicMessageSyntax ( )` `[inline]`

27.61.2.2 `virtual gdcmm::CryptographicMessageSyntax::~CryptographicMessageSyntax ( )` `[inline],[virtual]`

## 27.61.3 Member Function Documentation

27.61.3.1 `virtual bool gdcmm::CryptographicMessageSyntax::Decrypt ( char * output, size_t & outlen, const char * array, size_t len ) const` `[pure virtual]`

decrypt content from a CMS envelopedData structure

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), [gdcmm::CAPICryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.2 `virtual bool gdcmm::CryptographicMessageSyntax::Encrypt ( char * output, size_t & outlen, const char * array, size_t len ) const` `[pure virtual]`

create a CMS envelopedData structure

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), [gdcmm::CAPICryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.3 `virtual CipherTypes gdcmm::CryptographicMessageSyntax::GetCipherType ( ) const` `[pure virtual]`

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), [gdcmm::CAPICryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.4 `virtual bool gdcmm::CryptographicMessageSyntax::ParseCertificateFile ( const char * filename )` `[pure virtual]`

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), [gdcmm::CAPICryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.5 `virtual bool gdcmm::CryptographicMessageSyntax::ParseKeyFile ( const char * filename )` `[pure virtual]`

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), [gdcmm::CAPICryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.6 `virtual void gdcmm::CryptographicMessageSyntax::SetCipherType ( CipherTypes type )` `[pure virtual]`

Implemented in [gdcmm::OpenSSLP7CryptographicMessageSyntax](#), and [gdcmm::OpenSSLCryptographicMessageSyntax](#).

27.61.3.7 `virtual bool gdcm::CryptographicMessageSyntax::SetPassword ( const char * pass, size_t passLen ) [pure virtual]`

Implemented in [gdcm::OpenSSLP7CryptographicMessageSyntax](#), [gdcm::CAPICryptographicMessageSyntax](#), and [gdcm::OpenSSLCryptographicMessageSyntax](#).

The documentation for this class was generated from the following file:


- [gdcmCryptographicMessageSyntax.h](#)

## 27.62 gdcm::CSAElement Class Reference

Class to represent a CSA [Element](#).

```
#include <gdcmCSAElement.h>
```

Collaboration diagram for gdcm::CSAElement:

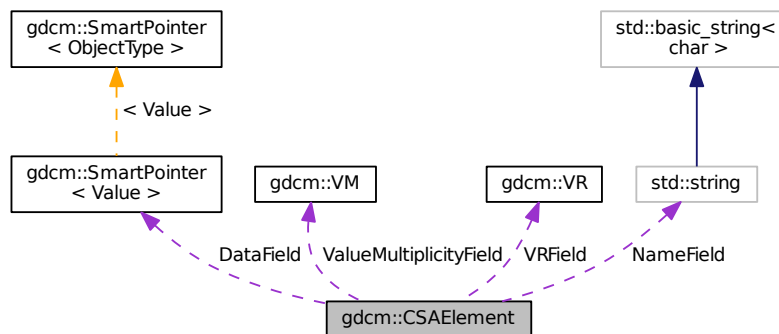

### Public Member Functions

- [CSAElement](#) (unsigned int kf=0)
- [CSAElement](#) (const [CSAElement](#) &\_val)
- const [ByteValue](#) \* [GetByteValue](#) () const
- unsigned int [GetKey](#) () const  
*Set/Get Key.*
- const char \* [GetName](#) () const  
*Set/Get Name.*
- unsigned int [GetNoOfItems](#) () const  
*Set/Get NoOfItems.*
- unsigned int [GetSyngoDT](#) () const  
*Set/Get SyngoDT.*
- [Value](#) const & [GetValue](#) () const  
*Set/Get Value (bytes array, SQ of items, SQ of fragments):*
- [Value](#) & [GetValue](#) ()

- const [VM](#) & [GetVM](#) () const  
*Set/Get VM.*
- [VR](#) const & [GetVR](#) () const  
*Set/Get VR.*
- bool [IsEmpty](#) () const  
*Check if CSA [Element](#) is empty.*
- bool [operator<](#) (const [CSAElement](#) &de) const
- [CSAElement](#) & [operator=](#) (const [CSAElement](#) &de)
- bool [operator==](#) (const [CSAElement](#) &de) const
- void [SetByteValue](#) (const char \*array, [VL](#) length)  
*Set.*
- void [SetKey](#) (unsigned int key)
- void [SetName](#) (const char \*name)
- void [SetNoOfItems](#) (unsigned int items)
- void [SetSyngoDT](#) (unsigned int syngodt)
- void [SetValue](#) ([Value](#) const &vl)
- void [SetVM](#) (const [VM](#) &vm)
- void [SetVR](#) ([VR](#) const &vr)

## Protected Types

- typedef [SmartPointer](#)< [Value](#) > [DataPtr](#)

## Protected Attributes

- [DataPtr](#) [DataField](#)
- unsigned int [KeyField](#)
- std::string [NameField](#)
- unsigned int [NoOfItemsField](#)
- unsigned int [SyngoDTField](#)
- [VM](#) [ValueMultiplicityField](#)
- [VR](#) [VRField](#)

## Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [CSAElement](#) &val)

### 27.62.1 Detailed Description

Class to represent a CSA [Element](#).

See also

[CSAHeader](#)

Examples:

[csa2img.cxx](#), and [MrProtocol.cxx](#).

## 27.62.2 Member Typedef Documentation

27.62.2.1 `typedef SmartPointer<Value> gdcm::CSAElement::DataPtr` `[protected]`

## 27.62.3 Constructor & Destructor Documentation

27.62.3.1 `gdcm::CSAElement::CSAElement ( unsigned int kf = 0 )` `[inline]`

27.62.3.2 `gdcm::CSAElement::CSAElement ( const CSAElement &_val )` `[inline]`

## 27.62.4 Member Function Documentation

27.62.4.1 `const ByteValue* gdcm::CSAElement::GetByteValue ( ) const` `[inline]`

Return the [Value](#) of [CSAElement](#) as a [ByteValue](#) (if possible)

Warning

: You need to check for NULL return value

Examples:

[MrProtocol.cxx](#).

27.62.4.2 `unsigned int gdcm::CSAElement::GetKey ( ) const` `[inline]`

Set/Get Key.

Referenced by operator<().

27.62.4.3 `const char* gdcm::CSAElement::GetName ( ) const` `[inline]`

Set/Get Name.

27.62.4.4 `unsigned int gdcm::CSAElement::GetNoOfItems ( ) const` `[inline]`

Set/Get NoOfItems.

27.62.4.5 `unsigned int gdcm::CSAElement::GetSyngoDT ( ) const` `[inline]`

Set/Get SyngoDT.

27.62.4.6 `Value const& gdcm::CSAElement::GetValue ( ) const` `[inline]`

Set/Get [Value](#) (bytes array, SQ of items, SQ of fragments):

Examples:

[csa2img.cxx](#).

27.62.4.7 **Value&** gdcm::CSAElement::GetValue ( ) [inline]

27.62.4.8 **const VM&** gdcm::CSAElement::GetVM ( ) **const** [inline]

Set/Get [VM](#).

27.62.4.9 **VR const&** gdcm::CSAElement::GetVR ( ) **const** [inline]

Set/Get [VR](#).

27.62.4.10 **bool** gdcm::CSAElement::IsEmpty ( ) **const** [inline]

Check if CSA [Element](#) is empty.

Examples:

[csa2img.cxx](#).

27.62.4.11 **bool** gdcm::CSAElement::operator< ( **const CSAElement & de** ) **const** [inline]

References [GetKey\(\)](#).

27.62.4.12 **CSAElement&** gdcm::CSAElement::operator= ( **const CSAElement & de** ) [inline]

References [DataField](#), [KeyField](#), [NameField](#), [NoOfItemsField](#), [SyngoDTField](#), [ValueMultiplicityField](#), and [VRField](#).

27.62.4.13 **bool** gdcm::CSAElement::operator== ( **const CSAElement & de** ) **const** [inline]

References [KeyField](#), [NameField](#), [SyngoDTField](#), [ValueMultiplicityField](#), and [VRField](#).

27.62.4.14 **void** gdcm::CSAElement::SetByteValue ( **const char \* array**, **VL length** ) [inline]

Set.

27.62.4.15 **void** gdcm::CSAElement::SetKey ( **unsigned int key** ) [inline]

27.62.4.16 **void** gdcm::CSAElement::SetName ( **const char \* name** ) [inline]

27.62.4.17 **void** gdcm::CSAElement::SetNoOfItems ( **unsigned int items** ) [inline]

27.62.4.18 **void** gdcm::CSAElement::SetSyngoDT ( **unsigned int syngodt** ) [inline]

27.62.4.19 **void** gdcm::CSAElement::SetValue ( **Value const & vl** ) [inline]

27.62.4.20 **void** gdcm::CSAElement::SetVM ( **const VM & vm** ) [inline]

27.62.4.21 **void** gdcm::CSAElement::SetVR ( **VR const & vr** ) [inline]

### 27.62.5 Friends And Related Function Documentation

27.62.5.1 `std::ostream& operator<< ( std::ostream & os, const CSAElement & val )` `[friend]`

### 27.62.6 Member Data Documentation

27.62.6.1 `DataPtr gdcM::CSAElement::DataField` `[protected]`

Referenced by `gdcM::operator<<()`, and `operator=()`.

27.62.6.2 `unsigned int gdcM::CSAElement::KeyField` `[protected]`

Referenced by `gdcM::operator<<()`, `operator=()`, and `operator==()`.

27.62.6.3 `std::string gdcM::CSAElement::NameField` `[protected]`

Referenced by `gdcM::operator<<()`, `operator=()`, and `operator==()`.

27.62.6.4 `unsigned int gdcM::CSAElement::NoOfItemsField` `[protected]`

Referenced by `gdcM::operator<<()`, and `operator=()`.

27.62.6.5 `unsigned int gdcM::CSAElement::SyngoDTField` `[protected]`

Referenced by `gdcM::operator<<()`, `operator=()`, and `operator==()`.

27.62.6.6 `VM gdcM::CSAElement::ValueMultiplicityField` `[protected]`

Referenced by `gdcM::operator<<()`, `operator=()`, and `operator==()`.

27.62.6.7 `VR gdcM::CSAElement::VRField` `[protected]`

Referenced by `gdcM::operator<<()`, `operator=()`, and `operator==()`.

The documentation for this class was generated from the following file:

- [gdcMCSAElement.h](#)

## 27.63 gdcM::CSAHeader Class Reference

Class for [CSAHeader](#).

```
#include <gdcMCSAHeader.h>
```



## Public Types

- enum [CSAHeaderType](#) {  
[UNKNOWN](#) = 0,  
[SV10](#),  
[NOMAGIC](#),  
[DATASET\\_FORMAT](#),  
[INTERFILE](#),  
[ZEROED\\_OUT](#) }

*Divers format of [CSAHeader](#) as found 'in the wild'.*

## Public Member Functions

- [CSAHeader](#) ()
- [~CSAHeader](#) ()
- bool [FindCSAElementByName](#) (const char \*name)
- const [CSAElement](#) & [GetCSAElementByName](#) (const char \*name)
- const [DataSet](#) & [GetDataSet](#) () const  
*Return the [DataSet](#) output (use only if Format == DATASET\_FORMAT)*
- [CSAHeaderType](#) [GetFormat](#) () const
- const char \* [GetInterfile](#) () const  
*Return the string output (use only if Format == Interfile)*
- bool [LoadFromDataElement](#) ([DataElement](#) const &de)  
*Decode the [CSAHeader](#) from element 'de'.*
- void [Print](#) (std::ostream &os) const  
*Print the [CSAHeader](#) (use only if Format == SV10 or NOMAGIC)*
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TSwap >  
const std::ostream & [Write](#) (std::ostream &os) const

## Static Public Member Functions

- static const [PrivateTag](#) & [GetCSADataInfo](#) ()
- static const [PrivateTag](#) & [GetCSAImageHeaderInfoTag](#) ()
- static const [PrivateTag](#) & [GetCSASeriesHeaderInfoTag](#) ()

## Protected Member Functions

- const [CSAElement](#) & [GetCSAEEnd](#) () const

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [CSAHeader](#) &d)

### 27.63.1 Detailed Description

Class for [CSAHeader](#).

SIEMENS store private information in tag (0x0029,0x10,"SIEMENS CSA HEADER") this class is meant for user wishing to access values stored within this private attribute. There are basically two main 'format' for this attribute : SV10/NOMAGIC and DATASET\_FORMAT SV10 and NOMAGIC are from a user prospective identical, see CSAHeader.xml for possible name / value stored in this format. DATASET\_FORMAT is in fact simply just another DICOM dataset (implicit) with -currently unknown- value. This can be only be printed for now.

#### Warning

Everything you do with this code is at your own risk, since decoding process was not written from specification documents.  
the API of this class might change.

**Todo** MrEvaProtocol in 29,1020 contains ^M that would be nice to get rid of on UNIX system...

See also

[PDBHeader](#)

External references: 5.1.3.2.4.1 MEDCOM History Information and 5.1.4.3 CSA Non-Image [Module](#) in [http://tamsinfo.toshiba.com/docrequest/pdf/E.Soft\\_v2.0.pdf](http://tamsinfo.toshiba.com/docrequest/pdf/E.Soft_v2.0.pdf)

Examples:

[csa2img.cxx](#), and [MrProtocol.cxx](#).

### 27.63.2 Member Enumeration Documentation

#### 27.63.2.1 enum gdcm::CSAHeader::CSAHeaderType

Divers format of [CSAHeader](#) as found 'in the wild'.

Enumerator

**UNKNOWN**  
**SV10**  
**NOMAGIC**  
**DATASET\_FORMAT**  
**INTERFILE**  
**ZEROED\_OUT**

### 27.63.3 Constructor & Destructor Documentation

27.63.3.1 `gdcm::CSAHeader::CSAHeader ( )` [`inline`]

27.63.3.2 `gdcm::CSAHeader::~~CSAHeader ( )` [`inline`]

### 27.63.4 Member Function Documentation

#### 27.63.4.1 bool gdcm::CSAHeader::FindCSAElementByName ( const char \* *name* )

Return true if the CSA element matching 'name' is found or not

Warning

Case Sensitive

Examples:

[csa2img.cxx](#), and [MrProtocol.cxx](#).

#### 27.63.4.2 static const PrivateTag& gdcm::CSAHeader::GetCSADataInfo ( ) [static]

Return the private tag used by SIEMENS to store the CSA Data Info This is: [PrivateTag](#)(0x0029,0x0010,"SIEMENS CSA NON-IMAGE");

#### 27.63.4.3 const CSAElement& gdcm::CSAHeader::GetCSAEnd ( ) const [protected]

#### 27.63.4.4 const CSAElement& gdcm::CSAHeader::GetCSAElementByName ( const char \* *name* )

Return the [CSAElement](#) corresponding to name 'name'

Warning

Case Sensitive

Examples:

[csa2img.cxx](#), and [MrProtocol.cxx](#).

#### 27.63.4.5 static const PrivateTag& gdcm::CSAHeader::GetCSAImageHeaderInfoTag ( ) [static]

Return the private tag used by SIEMENS to store the CSA [Image](#) Header This is: [PrivateTag](#)(0x0029,0x0010,"SIEMENS CSA HEADER");

Examples:

[csa2img.cxx](#), and [PublicDict.cxx](#).

#### 27.63.4.6 static const PrivateTag& gdcm::CSAHeader::GetCSASeriesHeaderInfoTag ( ) [static]

Return the private tag used by SIEMENS to store the CSA [Series](#) Header This is: [PrivateTag](#)(0x0029,0x0020,"SIEMENS CSA HEADER");

Examples:

[MrProtocol.cxx](#).

27.63.4.7 `const DataSet& gdcm::CSAHeader::GetDataSet ( ) const` `[inline]`

Return the [DataSet](#) output (use only if Format == DATASET\_FORMAT )

27.63.4.8 `CSAHeaderType gdcm::CSAHeader::GetFormat ( ) const`

return the format of the [CSAHeader](#) SV10 and NOMAGIC are equivalent.

27.63.4.9 `const char* gdcm::CSAHeader::GetInterfile ( ) const` `[inline]`

Return the string output (use only if Format == Interfile)

27.63.4.10 `bool gdcm::CSAHeader::LoadFromDataElement ( DataElement const & de )`

Decode the [CSAHeader](#) from element 'de'.

Examples:

[csa2img.cxx](#), and [MrProtocol.cxx](#).

27.63.4.11 `void gdcm::CSAHeader::Print ( std::ostream & os ) const`

Print the [CSAHeader](#) (use only if Format == SV10 or NOMAGIC)

Examples:

[csa2img.cxx](#).

Referenced by `gdcm::operator<<()`.

27.63.4.12 `template<typename TSwap > std::istream& gdcm::CSAHeader::Read ( std::istream & is )`

27.63.4.13 `template<typename TSwap > const std::ostream& gdcm::CSAHeader::Write ( std::ostream & os ) const`

## 27.63.5 Friends And Related Function Documentation

27.63.5.1 `std::ostream& operator<< ( std::ostream & _os, const CSAHeader & d )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmCSAHeader.h](#)

## 27.64 gdcm::CSAHeaderDict Class Reference

Class to represent a map of [CSAHeaderDictEntry](#).

```
#include <gdcmCSAHeaderDict.h>
```

## Public Types

- typedef MapCSAHeaderDictEntry::const\_iterator [ConstIterator](#)
- typedef MapCSAHeaderDictEntry::iterator [Iterator](#)
- typedef std::set< [CSAHeaderDictEntry](#) > [MapCSAHeaderDictEntry](#)

## Public Member Functions

- [CSAHeaderDict](#) ()
- void [AddCSAHeaderDictEntry](#) (const [CSAHeaderDictEntry](#) &de)
- [ConstIterator](#) [Begin](#) () const
- [ConstIterator](#) [End](#) () const
- const [CSAHeaderDictEntry](#) & [GetCSAHeaderDictEntry](#) (const char \*name) const
- bool [IsEmpty](#) () const

## Protected Member Functions

- void [LoadDefault](#) ()

## Friends

- class [Dicts](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [CSAHeaderDict](#) &\_val)

### 27.64.1 Detailed Description

Class to represent a map of [CSAHeaderDictEntry](#).

Examples:

[MrProtocol.cxx](#).

### 27.64.2 Member Typedef Documentation

27.64.2.1 typedef MapCSAHeaderDictEntry::const\_iterator gdcm::CSAHeaderDict::ConstIterator

27.64.2.2 typedef MapCSAHeaderDictEntry::iterator gdcm::CSAHeaderDict::Iterator

27.64.2.3 typedef std::set<CSAHeaderDictEntry> gdcm::CSAHeaderDict::MapCSAHeaderDictEntry

### 27.64.3 Constructor & Destructor Documentation

27.64.3.1 gdcm::CSAHeaderDict::CSAHeaderDict ( ) `[inline]`

### 27.64.4 Member Function Documentation

27.64.4.1 void gdcm::CSAHeaderDict::AddCSAHeaderDictEntry ( const CSAHeaderDictEntry & de ) `[inline]`

27.64.4.2 `ConstIterator gdcM::CSAHeaderDict::Begin ( ) const` `[inline]`

27.64.4.3 `ConstIterator gdcM::CSAHeaderDict::End ( ) const` `[inline]`

27.64.4.4 `const CSAHeaderDictEntry& gdcM::CSAHeaderDict::GetCSAHeaderDictEntry ( const char * name ) const`  
`[inline]`

Examples:

[MrProtocol.cxx](#).

27.64.4.5 `bool gdcM::CSAHeaderDict::IsEmpty ( ) const` `[inline]`

27.64.4.6 `void gdcM::CSAHeaderDict::LoadDefault ( )` `[protected]`

## 27.64.5 Friends And Related Function Documentation

27.64.5.1 `friend class Dicts` `[friend]`

27.64.5.2 `std::ostream& operator<< ( std::ostream & _os, const CSAHeaderDict & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcMCSAHeaderDict.h](#)

## 27.65 gdcM::CSAHeaderDictEntry Class Reference

Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcM::Tag](#) to the needed information.

```
#include <gdcMCSAHeaderDictEntry.h>
```

### Public Member Functions

- [CSAHeaderDictEntry](#) (const char \*name="", [VR](#) const &vr=[VR::INVALID](#), [VM](#) const &vm=[VM::VM0](#), const char \*desc="")
- const char \* [GetDescription](#) () const  
*Set/Get Description.*
- const char \* [GetName](#) () const  
*Set/Get Name.*
- const [VM](#) & [GetVM](#) () const  
*Set/Get VM.*
- const [VR](#) & [GetVR](#) () const  
*Set/Get VR.*
- bool [operator<](#) (const [CSAHeaderDictEntry](#) &entry) const
- void [SetDescription](#) (const char \*desc)
- void [SetName](#) (const char \*name)
- void [SetVM](#) ([VM](#) const &vm)
- void [SetVR](#) (const [VR](#) &vr)

## Friends

- `std::ostream & operator<< (std::ostream &_os, const CSAHeaderDictEntry &_val)`

### 27.65.1 Detailed Description

Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcm::Tag](#) to the needed information.

#### Note

bla TODO FIXME: Need a PublicCSAHeaderDictEntry...indeed [CSAHeaderDictEntry](#) has a notion of retired which does not exist in PrivateCSAHeaderDictEntry...

#### See also

[gdcm::Dict](#)

#### Examples:

[MrProtocol.cxx](#).

### 27.65.2 Constructor & Destructor Documentation

27.65.2.1 `gdcm::CSAHeaderDictEntry::CSAHeaderDictEntry ( const char * name = " ", VR const & vr = VR::INVALID, VM const & vm = VM::VMO, const char * desc = " " ) [inline]`

### 27.65.3 Member Function Documentation

27.65.3.1 `const char* gdcm::CSAHeaderDictEntry::GetDescription ( ) const [inline]`

Set/Get Description.

27.65.3.2 `const char* gdcm::CSAHeaderDictEntry::GetName ( ) const [inline]`

Set/Get Name.

Referenced by `operator<()`.

27.65.3.3 `const VM& gdcm::CSAHeaderDictEntry::GetVM ( ) const [inline]`

Set/Get [VM](#).

27.65.3.4 `const VR& gdcm::CSAHeaderDictEntry::GetVR ( ) const [inline]`

Set/Get [VR](#).

27.65.3.5 `bool gdcm::CSAHeaderDictEntry::operator< ( const CSAHeaderDictEntry & entry ) const [inline]`

References `GetName()`.

27.65.3.6 void gdcM::CSAHeaderDictEntry::SetDescription ( const char \* *desc* ) [inline]

27.65.3.7 void gdcM::CSAHeaderDictEntry::SetName ( const char \* *name* ) [inline]

27.65.3.8 void gdcM::CSAHeaderDictEntry::SetVM ( VM const & *vm* ) [inline]

27.65.3.9 void gdcM::CSAHeaderDictEntry::SetVR ( const VR & *vr* ) [inline]

## 27.65.4 Friends And Related Function Documentation

27.65.4.1 std::ostream& operator<< ( std::ostream & *\_os*, const CSAHeaderDictEntry & *\_val* ) [friend]


The documentation for this class was generated from the following file:

- [gdcMCSAHeaderDictEntry.h](#)

## 27.66 gdcM::CSAHeaderDictException Class Reference


```
#include <gdcMCSAHeaderDict.h>
```

Inheritance diagram for gdcM::CSAHeaderDictException:

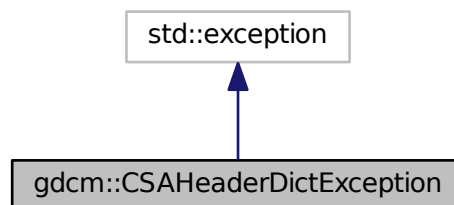



Collaboration diagram for gdcm::CSAHeaderDictException:

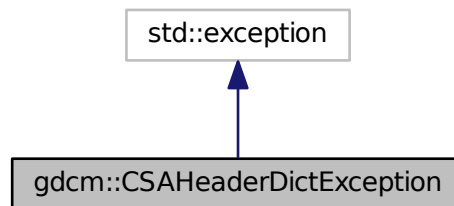

The documentation for this class was generated from the following file:


- [gdcmCSAHeaderDict.h](#)

## 27.67 gdcm::network::CStoreRQ Class Reference


[CStoreRQ](#) this file defines the messages for the cecho action.

```
#include <gdcmCStoreMessages.h>
```

Inheritance diagram for gdcm::network::CStoreRQ:

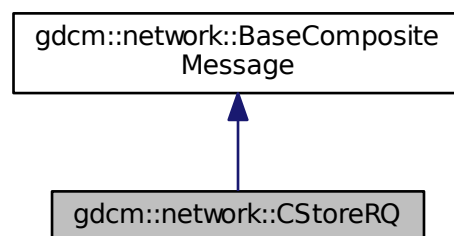

Collaboration diagram for `gdcm::network::CStoreRQ`:

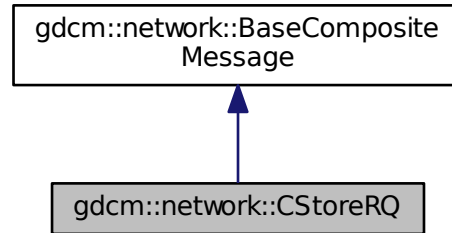

## Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDV` (const [ULConnection](#) &inConnection, const [File](#) &file)

### 27.67.1 Detailed Description

[CStoreRQ](#) this file defines the messages for the cecho action.

### 27.67.2 Member Function Documentation

**27.67.2.1** `std::vector<PresentationDataValue> gdcm::network::CStoreRQ::ConstructPDV ( const ULConnection &inConnection, const File &file )`

The documentation for this class was generated from the following file:


- [gdcmCStoreMessages.h](#)

## 27.68 `gdcm::network::CStoreRSP` Class Reference


[CStoreRSP](#) this file defines the messages for the cecho action.

```
#include <gdcmCStoreMessages.h>
```

Inheritance diagram for gdcmm::network::CStoreRSP:

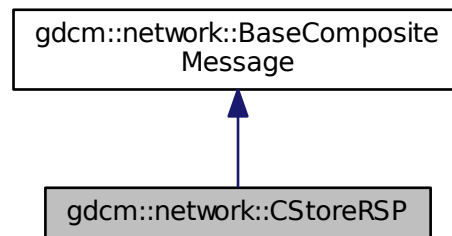

Collaboration diagram for gdcmm::network::CStoreRSP:

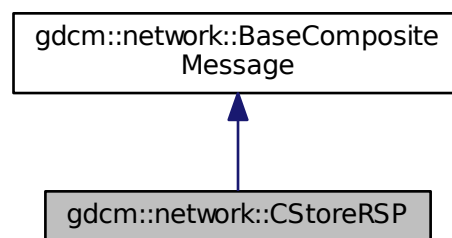

## Public Member Functions

- `std::vector< PresentationDataValue > ConstructPDV (const DataSet *inDataSet, const BasePDU *inPC)`

## 27.68.1 Detailed Description

[CStoreRSP](#) this file defines the messages for the cecho action.

## 27.68.2 Member Function Documentation

**27.68.2.1** `std::vector<PresentationDataValue> gdcmm::network::CStoreRSP::ConstructPDV ( const DataSet * inDataSet, const BasePDU * inPC )`

The documentation for this class was generated from the following file:


- [gdcmmCStoreMessages.h](#)

## 27.69 gdcM::Curve Class Reference


[Curve](#) class to handle element 50xx,3000 [Curve](#) Data WARNING: This is deprecated and lastly defined in PS 3.3 - 2004.

```
#include <gdcMCurve.h>
```

Inheritance diagram for gdcM::Curve:

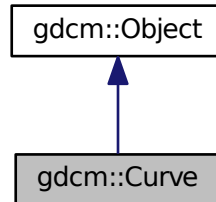

Collaboration diagram for gdcM::Curve:

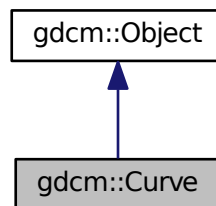

### Public Member Functions

- [Curve](#) ()
- [Curve](#) ([Curve](#) const &ov)
- [~Curve](#) ()
- void [Decode](#) (std::istream &is, std::ostream &os)
- void [GetAsPoints](#) (float \*array) const
- std::vector< unsigned short > const & [GetCurveDataDescriptor](#) () const
- unsigned short [GetDataValueRepresentation](#) () const
- unsigned short [GetDimensions](#) () const
- unsigned short [GetGroup](#) () const
- unsigned short [GetNumberOfPoints](#) () const

- const char \* [GetTypeOfData](#) () const
- const char \* [GetTypeOfDataDescription](#) () const
- bool [IsEmpty](#) () const
- void [Print](#) (std::ostream &) const
- void [SetCoordinateStartValue](#) (unsigned short v)
- void [SetCoordinateStepValue](#) (unsigned short v)
- void [SetCurve](#) (const char \*array, unsigned int length)
- void [SetCurveDataDescriptor](#) (const uint16\_t \*values, size\_t num)
- void [SetCurveDescription](#) (const char \*curvedescription)
- void [SetDataValueRepresentation](#) (unsigned short datavaluerepresentation)
- void [SetDimensions](#) (unsigned short dimensions)
- void [SetGroup](#) (unsigned short group)
- void [SetNumberOfPoints](#) (unsigned short numberofpoints)
- void [SetTypeOfData](#) (const char \*typeofdata)
- void [Update](#) (const [DataElement](#) &de)

### Static Public Member Functions

- static unsigned int [GetNumberOfCurves](#) ([DataSet](#) const &ds)

### Additional Inherited Members

#### 27.69.1 Detailed Description

[Curve](#) class to handle element 50xx,3000 [Curve](#) Data WARNING: This is deprecated and lastly defined in PS 3.3 - 2004.

Examples:

- GE\_DLX-8-MONO2-Multiframe-Jpeg\_Lossless.dcm
- GE\_DLX-8-MONO2-Multiframe.dcm
- gdcmSampleData/Philips\_Medical\_Images/integriss\_HV\_5000/xa\_integriss.dcm
- TOSHIBA-CurveData[1-3].dcm

#### 27.69.2 Constructor & Destructor Documentation

27.69.2.1 `gdcm::Curve::Curve ( )`

27.69.2.2 `gdcm::Curve::~~Curve ( )`

27.69.2.3 `gdcm::Curve::Curve ( Curve const & ov )`

#### 27.69.3 Member Function Documentation

27.69.3.1 `void gdcm::Curve::Decode ( std::istream & is, std::ostream & os )`

27.69.3.2 `void gdcm::Curve::GetAsPoints ( float * array ) const`

27.69.3.3 `std::vector<unsigned short> const& gdcM::Curve::GetCurveDataDescriptor ( ) const`

27.69.3.4 `unsigned short gdcM::Curve::GetDataValueRepresentation ( ) const`

27.69.3.5 `unsigned short gdcM::Curve::GetDimensions ( ) const`

27.69.3.6 `unsigned short gdcM::Curve::GetGroup ( ) const`

27.69.3.7 `static unsigned int gdcM::Curve::GetNumberOfCurves ( DataSet const & ds ) [static]`

27.69.3.8 `unsigned short gdcM::Curve::GetNumberOfPoints ( ) const`

27.69.3.9 `const char* gdcM::Curve::GetTypeOfData ( ) const`

27.69.3.10 `const char* gdcM::Curve::GetTypeOfDataDescription ( ) const`

27.69.3.11 `bool gdcM::Curve::IsEmpty ( ) const`

27.69.3.12 `void gdcM::Curve::Print ( std::ostream & ) const [virtual]`

Reimplemented from [gdcM::Object](#).

27.69.3.13 `void gdcM::Curve::SetCoordinateStartValue ( unsigned short v )`

27.69.3.14 `void gdcM::Curve::SetCoordinateStepValue ( unsigned short v )`

27.69.3.15 `void gdcM::Curve::SetCurve ( const char * array, unsigned int length )`

27.69.3.16 `void gdcM::Curve::SetCurveDataDescriptor ( const uint16_t * values, size_t num )`

27.69.3.17 `void gdcM::Curve::SetCurveDescription ( const char * curvedescription )`

27.69.3.18 `void gdcM::Curve::SetDataValueRepresentation ( unsigned short datavaluerepresentation )`

27.69.3.19 `void gdcM::Curve::SetDimensions ( unsigned short dimensions )`

27.69.3.20 `void gdcM::Curve::SetGroup ( unsigned short group )`

27.69.3.21 `void gdcM::Curve::SetNumberOfPoints ( unsigned short numberofpoints )`

27.69.3.22 `void gdcM::Curve::SetTypeOfData ( const char * typeofdata )`

27.69.3.23 `void gdcM::Curve::Update ( const DataElement & de )`

The documentation for this class was generated from the following file:


- [gdcMCurve.h](#)

## 27.70 gdcM::DataElement Class Reference


Class to represent a Data [Element](#) either Implicit or Explicit.

```
#include <gdcmDataElement.h>
```

Inheritance diagram for gdcm::DataElement:

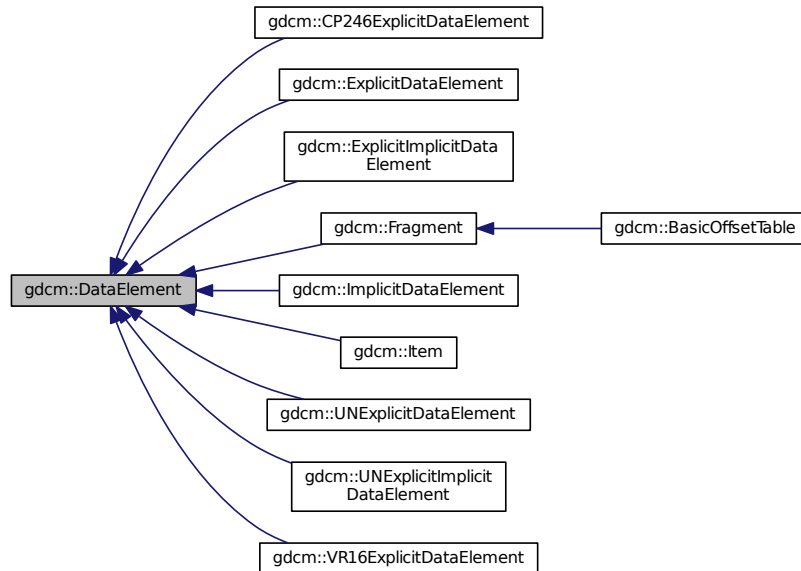

Collaboration diagram for gdcm::DataElement:

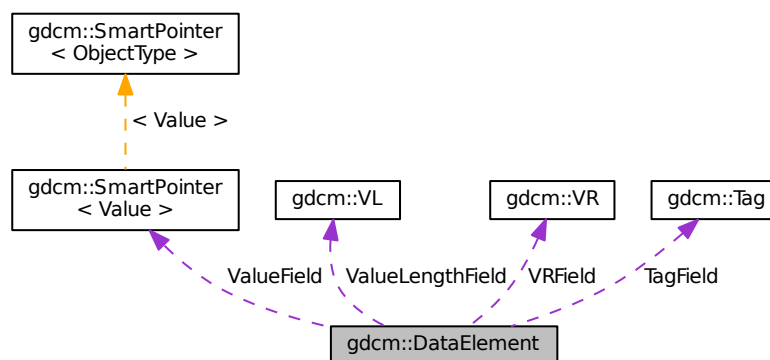

## Public Member Functions

- `DataElement` (const `Tag` &t=`Tag`(0), const `VL` &vl=0, const `VR` &vr=`VR::INVALID`)
- `DataElement` (const `DataElement` &\_val)
- void `Clear` ()

- Clear Data *Element* (make *Value* empty and invalidate *Tag* & *VR*)
- void **Empty** ()
  - Make Data *Element* empty (no *Value*)
- const **ByteValue** \* **GetByteValue** () const
- template<typename TDE >
  - **VL** **GetLength** () const
- const **SequenceOfFragments** \* **GetSequenceOfFragments** () const
- **SequenceOfFragments** \* **GetSequenceOfFragments** ()
- const **Tag** & **GetTag** () const
  - Get *Tag*.
- **Tag** & **GetTag** ()
- **Value** const & **GetValue** () const
  - Set/Get *Value* (bytes array, SQ of items, SQ of fragments):
- **Value** & **GetValue** ()
- **SmartPointer**< **SequenceOfItems** > **GetValueAsSQ** () const
- const **VL** & **GetVL** () const
  - Get *VL*.
- **VL** & **GetVL** ()
- **VR** const & **GetVR** () const
- bool **IsEmpty** () const
  - Check if Data *Element* is empty.
- bool **IsUndefinedLength** () const
  - return if *Value* Length if of undefined length
- bool **operator**< (const **DataElement** &de) const
- **DataElement** & **operator**= (const **DataElement** &de)
- bool **operator**== (const **DataElement** &de) const
- template<typename TDE , typename TSwap >
  - std::istream & **Read** (std::istream &is)
- template<typename TDE , typename TSwap >
  - std::istream & **ReadOrSkip** (std::istream &is, std::set< **Tag** > const &skiptags)
- template<typename TDE , typename TSwap >
  - std::istream & **ReadPreValue** (std::istream &is, std::set< **Tag** > const &skiptags)
- template<typename TDE , typename TSwap >
  - std::istream & **ReadValue** (std::istream &is, std::set< **Tag** > const &skiptags)
- template<typename TDE , typename TSwap >
  - std::istream & **ReadValueWithLength** (std::istream &is, **VL** &length, std::set< **Tag** > const &skiptags)
- template<typename TDE , typename TSwap >
  - std::istream & **ReadWithLength** (std::istream &is, **VL** &length)
- void **SetByteValue** (const char \*array, **VL** length)
- void **SetTag** (const **Tag** &t)
- void **SetValue** (**Value** const &vl)
- void **SetVL** (const **VL** &vl)
- void **SetVLToUndefined** ()
- void **SetVR** (**VR** const &vr)
- template<typename TDE , typename TSwap >
  - const std::ostream & **Write** (std::ostream &os) const

## Protected Types

- typedef **SmartPointer**< **Value** > **ValuePtr**



## Protected Member Functions

- void [SetValueFieldLength](#) (VL vl, bool readvalues)

## Protected Attributes

- Tag TagField
- ValuePtr ValueField
- VL ValueLengthField
- VR VRField

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [DataElement](#) &\_val)

### 27.70.1 Detailed Description

Class to represent a Data [Element](#) either Implicit or Explicit.

DATA ELEMENT: A unit of information as defined by a single entry in the data dictionary. An encoded Information Object Definition (IOD) [Attribute](#) that is composed of, at a minimum, three fields: a Data [Element Tag](#), a [Value](#) Length, and a [Value](#) Field. For some specific Transfer Syntaxes, a Data [Element](#) also contains a [VR](#) Field where the [Value](#) Representation of that Data [Element](#) is specified explicitly.

Design:

- A [DataElement](#) in GDCM always store [VL](#) ([Value](#) Length) on a 32 bits integer even when [VL](#) is 16 bits
- A [DataElement](#) always store the [VR](#) even for Implicit TS, in which case [VR](#) is defaulted to [VR::INVALID](#)
- For [Item](#) start/end (See 0xffff tags), [Value](#) is NULL

See also

[ExplicitDataElement](#) [ImplicitDataElement](#)

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllIVR.cxx](#), [GenFakelIdentifyFile.cxx](#), [GenFakelImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [rle2img.cxx](#), and [StreamReaderTest.cxx](#).

### 27.70.2 Member Typedef Documentation

27.70.2.1 `typedef SmartPointer<Value> gdcm::DataElement::ValuePtr` [protected]

### 27.70.3 Constructor & Destructor Documentation

27.70.3.1 `gdcm::DataElement::DataElement ( const Tag & t = Tag ( 0 ), const VL & vl = 0, const VR & vr = VR::INVALID )` `[inline]`

27.70.3.2 `gdcm::DataElement::DataElement ( const DataElement & _val )` `[inline]`

### 27.70.4 Member Function Documentation

27.70.4.1 `void gdcm::DataElement::Clear ( )` `[inline]`

Clear Data [Element](#) (make [Value](#) empty and invalidate [Tag](#) & [VR](#))

References `gdcm::VR::INVALID`.

Referenced by `gdcm::Item::Clear()`.

27.70.4.2 `void gdcm::DataElement::Empty ( )` `[inline]`

Make Data [Element](#) empty (no [Value](#))

27.70.4.3 `const ByteValue* gdcm::DataElement::GetByteValue ( ) const` `[inline]`

Return the [Value](#) of [DataElement](#) as a [ByteValue](#) (if possible)

Warning

: You need to check for NULL return value

Examples:

[DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetSubSequenceData.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcm::operator<<()`, `gdcm::Element< VR::OB, VM::VM1_n >::SetFromDataElement()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, `gdcm::Element< TVR, VM::VM1_n >::SetFromDataElement()`, and `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`.

27.70.4.4 `template<typename TDE> VL gdcm::DataElement::GetLength ( ) const` `[inline]`

27.70.4.5 `const SequenceOfFragments* gdcm::DataElement::GetSequenceOfFragments ( ) const`

Return the [Value](#) of [DataElement](#) as a Sequence Of Fragments (if possible)

Warning

: You need to check for NULL return value

Examples:

[FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), and [GetJPEGSamplePrecision.cxx](#).

27.70.4.6 **SequenceOfFragments\*** gdcm::DataElement::GetSequenceOfFragments ( )

27.70.4.7 **const Tag&** gdcm::DataElement::GetTag ( ) **const** [inline]

Get [Tag](#).

Examples:

[DumpGEMSMovieGroup.cxx](#), [DuplicatePCDE.cxx](#), [pmsct\\_rgb1.cxx](#), and [rle2img.cxx](#).

Referenced by [gdcm::CommandDataSet::Insert\(\)](#), [gdcm::FileMetaInformation::Insert\(\)](#), [gdcm::DataSet::Insert\(\)](#), [operator<\(\)](#), [gdcm::SequenceOfItems::Read\(\)](#), [gdcm::SequenceOfFragments::ReadValue\(\)](#), [gdcm::CommandDataSet::Replace\(\)](#), [gdcm::FileMetaInformation::Replace\(\)](#), [gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement\(\)](#), [gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement\(\)](#), and [gdcm::Attribute< Group, Element, TVR, VM::VM1\\_n >::SetFromDataElement\(\)](#).

27.70.4.8 **Tag&** gdcm::DataElement::GetTag ( ) [inline]

27.70.4.9 **Value const&** gdcm::DataElement::GetValue ( ) **const** [inline]

Set/Get [Value](#) (bytes array, SQ of items, SQ of fragments):

Examples:

[ReadAndDumpDICOMDIR.cxx](#).

Referenced by [gdcm::DataSet::InsertDataElement\(\)](#), [gdcm::Element< VR::OB, VM::VM1\\_n >::SetFromDataElement\(\)](#), and [gdcm::Element< TVR, VM::VM1\\_n >::SetFromDataElement\(\)](#).

27.70.4.10 **Value&** gdcm::DataElement::GetValue ( ) [inline]

27.70.4.11 **SmartPointer<SequenceOfItems>** gdcm::DataElement::GetValueAsSQ ( ) **const**

Interpret the [Value](#) stored in the [DataElement](#). This is more robust (but also more expensive) to call this function rather than the simplest form: [GetSequenceOfItems\(\)](#) It also return NULL when the [Value](#) is NOT of type [SequenceOfItems](#)

Warning

in case [GetSequenceOfItems\(\)](#) succeed the function return this value, otherwise it creates a new [SequenceOfItems](#), you should handle that in your case, for instance: `SmartPointer<SequenceOfItems> sqi = de.GetValueAsSQ();`

Examples:

[ChangeSequenceUltrasound.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [ExtractEncryptedContent.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [LargeVRDSExplicit.cxx](#), and [ReadAndDumpDICOMDIR.cxx](#).

27.70.4.12 **const VL&** gdcm::DataElement::GetVL ( ) **const** [inline]

Get [VL](#).

Referenced by [gdcm::DataSet::InsertDataElement\(\)](#), [gdcm::SequenceOfItems::Read\(\)](#), and [gdcm::SequenceOfFragments::ReadValue\(\)](#).

27.70.4.13 `VL& gdcM::DataElement::GetVL ( ) [inline]`

27.70.4.14 `VR const& gdcM::DataElement::GetVR ( ) const [inline]`

Get [VR](#) do not set [VR::SQ](#) on bytevalue data element

Examples:

[DuplicatePCDE.cxx](#), and [GenFakeIdentifyFile.cxx](#).

Referenced by `gdcM::Element< VR::OB, VM::VM1_n >::GetAsDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::GetAsDataElement()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement()`, `gdcM::Element< TVR, VM::VM1_n >::GetAsDataElement()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`, `gdcM::Element< VR::OB, VM::VM1_n >::SetFromDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, `gdcM::Element< TVR, VM::VM1_n >::SetFromDataElement()`, and `gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`.

27.70.4.15 `bool gdcM::DataElement::IsEmpty ( ) const [inline]`

Check if Data [Element](#) is empty.

Examples:

[DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [ELSCINT1WaveToText.cxx](#), [FixJAIBugJPEGLS.cxx](#), [pmsct\\_rgb1.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcM::DataSet::InsertDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`, `gdcM::Attribute< Group, Element, TVR, TVM >::SetFromDataSet()`, `gdcM::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet()`, and `gdcM::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataSet()`.

27.70.4.16 `bool gdcM::DataElement::IsUndefinedLength ( ) const [inline]`

return if [Value](#) Length if of undefined length

27.70.4.17 `bool gdcM::DataElement::operator< ( const DataElement & de ) const [inline]`

References `GetTag()`.

27.70.4.18 `DataElement& gdcM::DataElement::operator= ( const DataElement & de ) [inline]`

References `TagField`, `ValueField`, `ValueLengthField`, and `VRField`.

27.70.4.19 `bool gdcM::DataElement::operator== ( const DataElement & de ) const [inline]`

References `TagField`, `ValueField`, `ValueLengthField`, and `VRField`.

27.70.4.20 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::Read ( std::istream & is )`  
`[inline]`

27.70.4.21 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::ReadOrSkip ( std::istream & is,`  
`std::set< Tag > const & skiptags ) [inline]`

27.70.4.22 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::ReadPreValue ( std::istream & is,`  
`std::set< Tag > const & skiptags ) [inline]`

27.70.4.23 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::ReadValue ( std::istream & is,`  
`std::set< Tag > const & skiptags ) [inline]`

27.70.4.24 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::ReadValueWithLength ( std::istream &`  
`is, VL & length, std::set< Tag > const & skiptags ) [inline]`

27.70.4.25 `template<typename TDE , typename TSwap > std::istream& gdcm::DataElement::ReadWithLength ( std::istream & is,`  
`VL & length ) [inline]`

27.70.4.26 `void gdcm::DataElement::SetByteValue ( const char * array, VL length ) [inline]`

Set the byte value

#### Warning

user need to read DICOM standard for an understanding of:

- even padding
- \0 vs space padding By default even padding is achieved using \0 regardless of the of [VR](#)

#### Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FileChangeTS.cs](#), [FileChangeTSLossy.cs](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenFakeImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [NewSequence.cs](#), and [StreamImageReaderTest.cxx](#).

Referenced by `gdcm::Element< VR::OB, VM::VM1_n >::GetAsDataElement()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`, `gdcm::Element< TVR, VM::VM1_n >::GetAsDataElement()`, and `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`.

27.70.4.27 `void gdcm::DataElement::SetTag ( const Tag & t ) [inline]`

Set [Tag](#) Use with cautious (need to match Part 6)

#### Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenFakeIdentifyFile.cxx](#), and [GetSubSequenceData.cxx](#).

27.70.4.28 `void gdcm::DataElement::SetValue ( Value const & v ) [inline]`

**Warning**

you need to set the ValueLengthField explicitly

**Examples:**

[DuplicatePCDE.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), and [GenSeqs.cxx](#).

References `gdcmm::Value::GetLength()`.

**27.70.4.29** `void gdcmm::DataElement::SetValueFieldLength ( VL vl, bool readvalues )` [protected]

**27.70.4.30** `void gdcmm::DataElement::SetVL ( const VL & vl )` [inline]

Set [VL](#) Use with cautious (need to match Part 6), advanced user only

**See also**

[SetByteValue](#)

**27.70.4.31** `void gdcmm::DataElement::SetVLToUndefined ( )`

**Examples:**

[Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllIVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), and [NewSequence.cs](#).

**27.70.4.32** `void gdcmm::DataElement::SetVR ( VR const & vr )` [inline]

Set [VR](#) Use with cautious (need to match Part 6), advanced user only

**Precondition**

vr is a [VR::VRALL](#) (not a dual one such as OB\_OW)

**Examples:**

[Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAI BugJPEGLS.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [NewSequence.cs](#), and [StreamImageReaderTest.cxx](#).

References `gdcmm::VR::IsVRFile()`.

Referenced by `gdcmm::Element< VR::OB, VM::VM1_n >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, T < VR, TVM >::GetAsDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::GetAsDataElement()`, `gdcmm::Element< TVR, VM::VM1_n >::GetAsDataElement()`, and `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::GetAsDataElement()`.

27.70.4.33 `template<typename TDE , typename TSwap > const std::ostream& gdcm::DataElement::Write ( std::ostream & os )`  
`const [inline]`

## 27.70.5 Friends And Related Function Documentation

27.70.5.1 `std::ostream& operator<< ( std::ostream & _os, const DataElement & _val ) [friend]`

## 27.70.6 Member Data Documentation

27.70.6.1 `Tag gdcm::DataElement::TagField [protected]`

Referenced by `gdcm::operator<<()`, `operator=()`, and `operator==()`.

27.70.6.2 `ValuePtr gdcm::DataElement::ValueField [protected]`

Referenced by `gdcm::operator<<()`, `operator=()`, and `operator==()`.

27.70.6.3 `VL gdcm::DataElement::ValueLengthField [protected]`

Referenced by `gdcm::operator<<()`, `operator=()`, and `operator==()`.

27.70.6.4 `VR gdcm::DataElement::VRField [protected]`

Referenced by `gdcm::operator<<()`, `operator=()`, and `operator==()`.


The documentation for this class was generated from the following file:

- [gdcmDataElement.h](#)


## 27.71 gdcm::DataElementException Class Reference

```
#include <gdcmDataSet.h>
```

Inheritance diagram for `gdcm::DataElementException`:

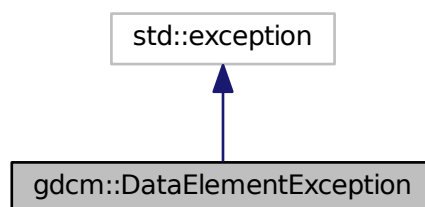

Collaboration diagram for `gdcm::DataElementException`:

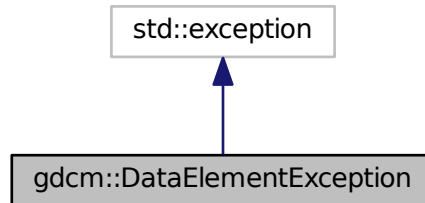

The documentation for this class was generated from the following file:


- [gdcmDataSet.h](#)

## 27.72 `gdcm::DataEvent` Class Reference

[DataEvent](#).


```
#include <gdcmDataEvent.h>
```

Inheritance diagram for `gdcm::DataEvent`:

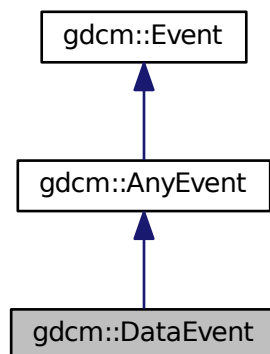



Collaboration diagram for gdcm::DataEvent:

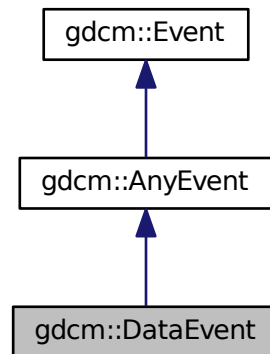

## Public Types

- typedef [DataEvent Self](#)
- typedef [AnyEvent Superclass](#)

## Public Member Functions

- [DataEvent](#) (const char \*bytes=0, size\_t len=0)
- [DataEvent](#) (const [Self](#) &s)
- virtual [~DataEvent](#) ()
- virtual bool [CheckEvent](#) (const ::gdcm::Event \*e) const
- const char \* [GetData](#) () const
- size\_t [GetDataLength](#) () const
- virtual const char \* [GetEventName](#) () const
- virtual ::gdcm::Event \* [MakeObject](#) () const
- void [SetData](#) (const char \*bytes, size\_t len)

### 27.72.1 Detailed Description

[DataEvent](#).

### 27.72.2 Member Typedef Documentation

27.72.2.1 typedef `DataEvent` `gdcm::DataEvent::Self`

27.72.2.2 typedef `AnyEvent` `gdcm::DataEvent::Superclass`

### 27.72.3 Constructor & Destructor Documentation

27.72.3.1 `gdcm::DataEvent::DataEvent ( const char * bytes = 0, size_t len = 0 ) [inline]`

27.72.3.2 `virtual gdcm::DataEvent::~~DataEvent ( ) [inline],[virtual]`

27.72.3.3 `gdcm::DataEvent::DataEvent ( const Self & s ) [inline]`

## 27.72.4 Member Function Documentation

27.72.4.1 `virtual bool gdcm::DataEvent::CheckEvent ( const ::gdcm::Event * e ) const [inline],[virtual]`

27.72.4.2 `const char* gdcm::DataEvent::GetData ( ) const [inline]`

27.72.4.3 `size_t gdcm::DataEvent::GetDataLength ( ) const [inline]`

27.72.4.4 `virtual const char* gdcm::DataEvent::GetEventName ( ) const [inline],[virtual]`

Return the StringName associated with the event.

Implements [gdcm::Event](#).

27.72.4.5 `virtual ::gdcm::Event* gdcm::DataEvent::MakeObject ( ) const [inline],[virtual]`

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implements [gdcm::Event](#).

27.72.4.6 `void gdcm::DataEvent::SetData ( const char * bytes, size_t len ) [inline]`

The documentation for this class was generated from the following file:


- [gdcmDataEvent.h](#)

## 27.73 gdcm::DataSet Class Reference

Class to represent a Data Set (which contains Data Elements) A Data Set represents an instance of a real world Information [Object](#).

```
#include <gdcmDataSet.h>
```

Inheritance diagram for gdcmm::DataSet:

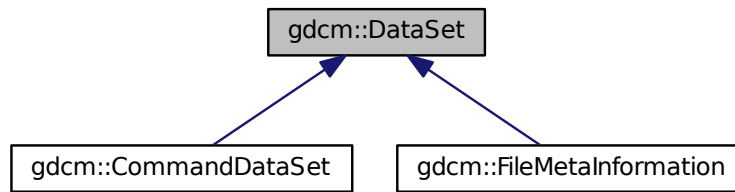

## Public Types

- typedef DataSet::const\_iterator [ConstIterator](#)
- typedef std::set< [DataElement](#) > [DataElementSet](#)
- typedef DataSet::iterator [Iterator](#)
- typedef DataSet::size\_type [SizeType](#)

## Public Member Functions

- [ConstIterator Begin](#) () const
- [Iterator Begin](#) ()
- void [Clear](#) ()
- template<typename TDE >  
unsigned int [ComputeGroupLength](#) ([Tag](#) const &tag) const
- [ConstIterator End](#) () const
- [Iterator End](#) ()
- bool [FindDataElement](#) (const [PrivateTag](#) &t) const  
*Look up if private tag 't' is present in the dataset:*
- bool [FindDataElement](#) (const [Tag](#) &t) const
- const [DataElement](#) & [FindNextDataElement](#) (const [Tag](#) &t) const
- const [DataElement](#) & [GetDataElement](#) (const [Tag](#) &t) const
- const [DataElement](#) & [GetDataElement](#) (const [PrivateTag](#) &t) const  
*Return the dataelement.*
- const [DataElementSet](#) & [GetDES](#) () const
- [DataElementSet](#) & [GetDES](#) ()
- template<typename TDE >  
[VL GetLength](#) () const
- [MediaStorage GetMediaStorage](#) () const
- std::string [GetPrivateCreator](#) (const [Tag](#) &t) const  
*Return the private creator of the private tag 't':*
- void [Insert](#) (const [DataElement](#) &de)
- bool [IsEmpty](#) () const  
*Returns if the dataset is empty.*
- const [DataElement](#) & [operator\(\)](#) (uint16\_t group, uint16\_t element) const

- [DataSet](#) & [operator=](#) ([DataSet](#) const &val)
- const [DataElement](#) & [operator\[\]](#) (const [Tag](#) &t) const
- void [Print](#) (std::ostream &os, std::string const &indent="") const
- template<typename TDE , typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TDE , typename TSwap >  
std::istream & [ReadNested](#) (std::istream &is)
- template<typename TDE , typename TSwap >  
std::istream & [ReadSelectedPrivateTags](#) (std::istream &is, const std::set< [PrivateTag](#) > &tags, bool readvalues=true)
- template<typename TDE , typename TSwap >  
std::istream & [ReadSelectedPrivateTagsWithLength](#) (std::istream &is, const std::set< [PrivateTag](#) > &tags, [VL](#) &length, bool readvalues=true)
- template<typename TDE , typename TSwap >  
std::istream & [ReadSelectedTags](#) (std::istream &is, const std::set< [Tag](#) > &tags, bool readvalues=true)
- template<typename TDE , typename TSwap >  
std::istream & [ReadSelectedTagsWithLength](#) (std::istream &is, const std::set< [Tag](#) > &tags, [VL](#) &length, bool readvalues=true)
- template<typename TDE , typename TSwap >  
std::istream & [ReadUpToTag](#) (std::istream &is, const [Tag](#) &t, std::set< [Tag](#) > const &skiptags)
- template<typename TDE , typename TSwap >  
std::istream & [ReadUpToTagWithLength](#) (std::istream &is, const [Tag](#) &t, std::set< [Tag](#) > const &skiptags, [VL](#) &length)
- template<typename TDE , typename TSwap >  
std::istream & [ReadWithLength](#) (std::istream &is, [VL](#) &length)
- [SizeType Remove](#) (const [Tag](#) &tag)  
*Completely remove a dataelement from the dataset.*
- void [Replace](#) (const [DataElement](#) &de)  
*Replace a dataelement with another one.*
- void [ReplaceEmpty](#) (const [DataElement](#) &de)  
*Only replace a DICOM attribute when it is missing or empty.*
- [SizeType Size](#) () const
- template<typename TDE , typename TSwap >  
std::ostream const & [Write](#) (std::ostream &os) const

## Protected Member Functions

- [Tag ComputeDataElement](#) (const [PrivateTag](#) &t) const
- const [DataElement](#) & [GetDEEnd](#) () const
- void [InsertDataElement](#) (const [DataElement](#) &de)

## Friends

- class [CSAHeader](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [DataSet](#) &val)

### 27.73.1 Detailed Description

Class to represent a Data Set (which contains Data Elements) A Data Set represents an instance of a real world Information [Object](#).

#### Note

DATA SET: Exchanged information consisting of a structured set of [Attribute](#) values directly or indirectly related to Information Objects. The value of each [Attribute](#) in a Data Set is expressed as a Data [Element](#). A collection of Data Elements ordered by increasing Data [Element Tag](#) number that is an encoding of the values of Attributes of a real world object.

Implementation note. If one do: [DataSet](#) ds; ds.SetLength(0); ds.Read(is); setting length to 0 actually means try to read is as if it was a root [DataSet](#). Other value are undefined (nested dataset with undefined length) or defined length (different from 0) means nested dataset with defined length.

#### Warning

a [DataSet](#) does not have a Transfer Syntax type, only a [File](#) does.

#### Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [gdcmrtpionplan.cxx](#), [gdcmrtpionplan.cxx](#), [GenAllIVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [rle2img.cxx](#), [SortImage.cxx](#), [StreamImageReaderTest.cxx](#), and [VolumeSorter.cxx](#).

### 27.73.2 Member Typedef Documentation

27.73.2.1 `typedef DataElementSet::const_iterator gdcm::DataSet::ConstIterator`

27.73.2.2 `typedef std::set<DataElement> gdcm::DataSet::DataElementSet`

27.73.2.3 `typedef DataElementSet::iterator gdcm::DataSet::Iterator`

27.73.2.4 `typedef DataElementSet::size_type gdcm::DataSet::SizeType`

### 27.73.3 Member Function Documentation

27.73.3.1 `ConstIterator gdcm::DataSet::Begin ( ) const` `[inline]`

#### Examples:

[DiffFile.cxx](#), [DumpGEMSMovieGroup.cxx](#), and [DuplicatePCDE.cxx](#).

27.73.3.2 `Iterator gdcm::DataSet::Begin ( )` `[inline]`

27.73.3.3 `void gdcm::DataSet::Clear ( )` `[inline]`

Referenced by `gdcm::Item::Read()`.

27.73.3.4 **Tag** `gdcm::DataSet::ComputeDataElement ( const PrivateTag & t ) const` `[protected]`

27.73.3.5 `template<typename TDE > unsigned int gdcm::DataSet::ComputeGroupLength ( Tag const & tag ) const` `[inline]`

References `gdcm::Tag::GetElement()`, and `gdcm::Tag::GetGroup()`.

27.73.3.6 **ConstIterator** `gdcm::DataSet::End ( ) const` `[inline]`

Examples:

[DiffFile.cxx](#), [DumpGEMSMovieGroup.cxx](#), and [DuplicatePCDE.cxx](#).

27.73.3.7 **Iterator** `gdcm::DataSet::End ( )` `[inline]`

27.73.3.8 **bool** `gdcm::DataSet::FindDataElement ( const PrivateTag & t ) const`

Look up if private tag 't' is present in the dataset:

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [csa2img.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [gdcmrptionplan.cxx](#), [gdcmrtpian.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [LargeVRDSExplicit.cxx](#), [MrProtocol.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataSet()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet()`, and `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataSet()`.

27.73.3.9 **bool** `gdcm::DataSet::FindDataElement ( const Tag & t ) const` `[inline]`

27.73.3.10 **const DataElement&** `gdcm::DataSet::FindNextDataElement ( const Tag & t ) const` `[inline]`

Examples:

[DuplicatePCDE.cxx](#).

27.73.3.11 **const DataElement&** `gdcm::DataSet::GetDataElement ( const Tag & t ) const` `[inline]`

Return the [DataElement](#) with [Tag](#) 't'

Warning

: This only search at the 'root level' of the [DataSet](#)

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [csa2img.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAI BugJPEGLS.cxx](#), [gdcmrptionplan.cxx](#), [gdcmrtpian.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

Referenced by `gdcm::Attribute< Group, Element, TVR, TVM >::Set()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::Set()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::Set()`, `gdcm::Attribute< Group, Element, TVR, TVM >::SetFromDataSet()`, `gdcm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataSet()`, and `gdcm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataSet()`.

**27.73.3.12** `const DataElement& gdcm::DataSet::GetDataElement ( const PrivateTag & t ) const`

Return the dataelement.

**27.73.3.13** `const DataElement& gdcm::DataSet::GetDEEnd ( ) const` `[protected]`

**27.73.3.14** `const DataElementSet& gdcm::DataSet::GetDES ( ) const` `[inline]`

Examples:

[ReadAndDumpDICOMDIR.cxx](#).

**27.73.3.15** `DataElementSet& gdcm::DataSet::GetDES ( )` `[inline]`

**27.73.3.16** `template<typename TDE > VL gdcm::DataSet::GetLength ( ) const` `[inline]`

**27.73.3.17** `MediaStorage gdcm::DataSet::GetMediaStorage ( ) const`

**27.73.3.18** `std::string gdcm::DataSet::GetPrivateCreator ( const Tag & t ) const`

Return the private creator of the private tag 't':

Examples:

[DuplicatePCDE.cxx](#).

**27.73.3.19** `void gdcm::DataSet::Insert ( const DataElement & de )` `[inline]`

Insert a [DataElement](#) in the [DataSet](#).

Warning

: [Tag](#) need to be  $\geq 0x8$  to be considered valid data element

Examples:

[CreateJIPIDataset.cxx](#), [DuplicatePCDE.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [NewSequence.cs](#), and [StreamImageReaderTest.cxx](#).

References `gdcmErrorMacro`, `gdcm::Tag::GetGroup()`, and `gdcm::DataElement::GetTag()`.

**27.73.3.20** `void gdcm::DataSet::InsertDataElement ( const DataElement & de )` `[inline]`, `[protected]`

References `gdcmWarningMacro`, `gdcm::Value::GetLength()`, `gdcm::DataElement::GetValue()`, `gdcm::DataElement::GetVL()`, and `gdcm::DataElement::IsEmpty()`.

27.73.3.21 `bool gdcmm::DataSet::IsEmpty ( ) const [inline]`

Returns if the dataset is empty.

Referenced by `gdcmm::Item::Read()`.

27.73.3.22 `const DataElement& gdcmm::DataSet::operator() ( uint16_t group, uint16_t element ) const [inline]`

27.73.3.23 `DataSet& gdcmm::DataSet::operator= ( DataSet const & val ) [inline]`

27.73.3.24 `const DataElement& gdcmm::DataSet::operator[] ( const Tag & t ) const [inline]`

27.73.3.25 `void gdcmm::DataSet::Print ( std::ostream & os, std::string const & indent = " " ) const [inline]`

Referenced by `gdcmm::operator<<()`.

27.73.3.26 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::Read ( std::istream & is )`

27.73.3.27 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadNested ( std::istream & is )`

27.73.3.28 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadSelectedPrivateTags ( std::istream & is, const std::set< PrivateTag > & tags, bool readvalues = true )`

27.73.3.29 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadSelectedPrivateTagsWithLength ( std::istream & is, const std::set< PrivateTag > & tags, VL & length, bool readvalues = true )`

27.73.3.30 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadSelectedTags ( std::istream & is, const std::set< Tag > & tags, bool readvalues = true )`

27.73.3.31 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadSelectedTagsWithLength ( std::istream & is, const std::set< Tag > & tags, VL & length, bool readvalues = true )`

27.73.3.32 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadUpToTag ( std::istream & is, const Tag & t, std::set< Tag > const & skiptags )`

27.73.3.33 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadUpToTagWithLength ( std::istream & is, const Tag & t, std::set< Tag > const & skiptags, VL & length )`

27.73.3.34 `template<typename TDE , typename TSwap > std::istream& gdcmm::DataSet::ReadWithLength ( std::istream & is, VL & length )`

27.73.3.35 `SizeType gdcmm::DataSet::Remove ( const Tag & tag ) [inline]`

Completely remove a dataelement from the dataset.

Examples:

[ClinicalTrialIdentificationWorkflow.cs](#), [GenFakeIdentifyFile.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [pmsct\\_rgb1.cxx](#), [ReformatFile.cs](#), [rle2img.cxx](#), and [StandardizeFiles.cs](#).



27.73.3.36 void gdcm::DataSet::Replace ( const DataElement & *de* ) [inline]

Replace a dataelement with another one.

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateFakeRTDOSE.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAI↵  
BugJPEGLS.cxx](#), [GenFakeIdentifyFile.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [Large↵  
VRDSExplicit.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), and [rle2img.cxx](#).

27.73.3.37 void gdcm::DataSet::ReplaceEmpty ( const DataElement & *de* ) [inline]

Only replace a DICOM attribute when it is missing or empty.

27.73.3.38 SizeType gdcm::DataSet::Size ( ) const [inline]

Examples:

[DumpGEMSMovieGroup.cxx](#).

Referenced by `gdcm::SequenceOfItems::Read()`.

27.73.3.39 template<typename TDE , typename TSwap > std::ostream const& gdcm::DataSet::Write ( std::ostream & *os* ) const

## 27.73.4 Friends And Related Function Documentation

27.73.4.1 friend class CSAHeader [friend]

27.73.4.2 std::ostream& operator<< ( std::ostream & *\_os*, const DataSet & *val* ) [friend]

The documentation for this class was generated from the following file:


- [gdcmDataSet.h](#)

## 27.74 gdcm::DataSetEvent Class Reference


**DataSetEvent** Special type of event triggered during the [DataSet](#) store/move process.

```
#include <gdcmDataSetEvent.h>
```

Inheritance diagram for `gdcm::DataSetEvent`:

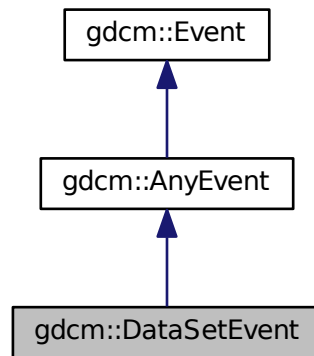

Collaboration diagram for `gdcm::DataSetEvent`:

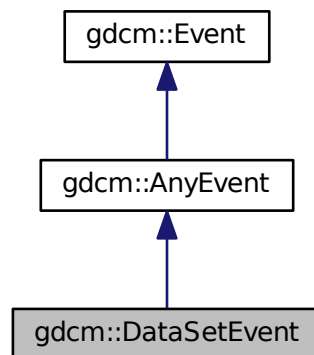

## Public Types

- typedef [DataSetEvent](#) `Self`
- typedef [AnyEvent](#) `Superclass`

## Public Member Functions

- [DataSetEvent](#) ([DataSet](#) const \*ds=NULL)
- [DataSetEvent](#) (const [Self](#) &s)

- virtual [~DataSetEvent](#) ()
- virtual bool [CheckEvent](#) (const ::gdcm::Event \*e) const
- [DataSet](#) const & [GetDataSet](#) () const
- virtual const char \* [GetEventName](#) () const
- virtual ::gdcm::Event \* [MakeObject](#) () const

### 27.74.1 Detailed Description

[DataSetEvent](#) Special type of event triggered during the [DataSet](#) store/move process.

See also

### 27.74.2 Member Typedef Documentation

27.74.2.1 typedef [DataSetEvent](#) gdcm::DataSetEvent::Self

27.74.2.2 typedef [AnyEvent](#) gdcm::DataSetEvent::Superclass

### 27.74.3 Constructor & Destructor Documentation

27.74.3.1 gdcm::DataSetEvent::DataSetEvent ( [DataSet](#) const \* *ds* = NULL ) [inline]

27.74.3.2 virtual gdcm::DataSetEvent::~~DataSetEvent ( ) [inline],[virtual]

27.74.3.3 gdcm::DataSetEvent::DataSetEvent ( const Self & *s* ) [inline]

### 27.74.4 Member Function Documentation

27.74.4.1 virtual bool gdcm::DataSetEvent::CheckEvent ( const ::gdcm::Event \* *e* ) const [inline],[virtual]

27.74.4.2 [DataSet](#) const& gdcm::DataSetEvent::GetDataSet ( ) const [inline]

27.74.4.3 virtual const char\* gdcm::DataSetEvent::GetEventName ( ) const [inline],[virtual]

Return the StringName associated with the event.

Implements [gdcm::Event](#).

27.74.4.4 virtual ::gdcm::Event\* gdcm::DataSetEvent::MakeObject ( ) const [inline],[virtual]

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implements [gdcm::Event](#).

The documentation for this class was generated from the following file:

- [gdcmDataSetEvent.h](#)

## 27.75 gdcm::DataSetHelper Class Reference

[DataSetHelper](#) (internal class, not intended for user level)

```
#include <gdcmDataSetHelper.h>
```

### Static Public Member Functions

- static [VR ComputeVR](#) ([File](#) const &file, [DataSet](#) const &ds, const [Tag](#) &tag)

### 27.75.1 Detailed Description

[DataSetHelper](#) (internal class, not intended for user level)

### 27.75.2 Member Function Documentation

27.75.2.1 static VR [gdcm::DataSetHelper::ComputeVR](#) ( [File](#) const & *file*, [DataSet](#) const & *ds*, const [Tag](#) & *tag* ) [static]

ds -> current dataset, which is not the same as the root dataset return [VR::INVALID](#) in case of error

The documentation for this class was generated from the following file:


- [gdcmDataSetHelper.h](#)

## 27.76 gdcm::Decoder Class Reference

[Decoder](#).

```
#include <gdcmDecoder.h>
```

Inheritance diagram for [gdcm::Decoder](#):

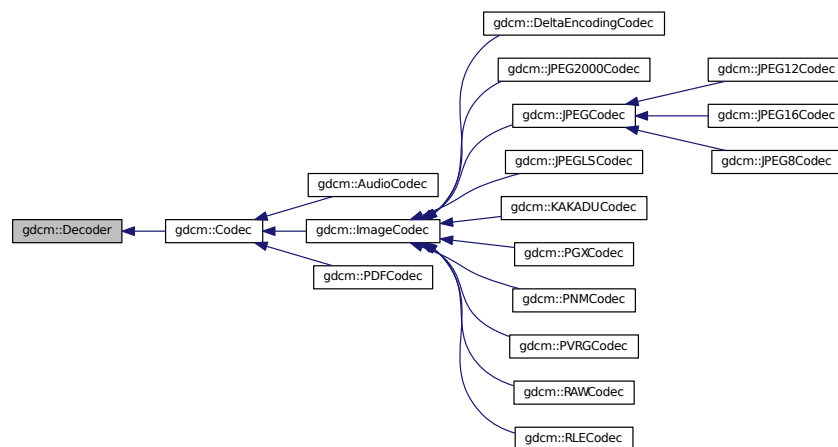

## Public Member Functions

- virtual [~Decoder](#) ()
- virtual bool [CanDecode](#) ([TransferSyntax](#) const &) const =0  
*Return whether this decoder support this transfer syntax (can decode it)*
- virtual bool [Decode](#) ([DataElement](#) const &, [DataElement](#) &)  
*Decode.*

## Protected Member Functions

- virtual bool [DecodeByStreams](#) (std::istream &, std::ostream &)

### 27.76.1 Detailed Description

[Decoder](#).

### 27.76.2 Constructor & Destructor Documentation

27.76.2.1 virtual gdcm::Decoder::~~Decoder ( ) [inline],[virtual]

### 27.76.3 Member Function Documentation

27.76.3.1 virtual bool gdcm::Decoder::CanDecode ( [TransferSyntax](#) const & ) const [pure virtual]

Return whether this decoder support this transfer syntax (can decode it)

Implemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::ImageCodec](#), [gdcm::JPEGLSCo](#)  
[Codec](#), [gdcm::PNMCodec](#), [gdcm::RAWCodec](#), [gdcm::AudioCodec](#), [gdcm::PDFCodec](#), [gdcm::PGXCodec](#), and [gdcm::KAKADUCodec](#).

27.76.3.2 virtual bool gdcm::Decoder::Decode ( [DataElement](#) const &, [DataElement](#) & ) [inline],[virtual]

Decode.

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::JPEGLSCo](#)  
[Codec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::ImageCodec](#), [gdcm::DeltaEncodingCodec](#), [gdcm::KAKADUCodec](#), [gdcm::RAWCodec](#), [gdcm::AudioCodec](#), and [gdcm::PDFCodec](#).

27.76.3.3 virtual bool gdcm::Decoder::DecodeByStreams ( std::istream &, std::ostream & ) [inline],[protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::RLECodec](#), [gdcm::ImageCodec](#), [gdcm::RAWCodec](#), [gdcm::JPEG12Codec](#), [gdcm::JPEG16Codec](#), and [gdcm::JPEG8Codec](#).

The documentation for this class was generated from the following file:

- [gdcmDecoder.h](#)

## 27.77 gdcm::DefinedTerms Class Reference

Defined Terms are used when the specified explicit Values may be extended by implementors to include additional new Values. These new Values shall be specified in the Conformance Statement (see PS 3.2) and shall not have the same meaning as currently defined Values in this standard. A Data [Element](#) with Defined Terms that does not contain a [Value](#) equivalent to one of the Values currently specified in this standard shall not be considered to have an invalid value. Note: Interpretation [Type](#) ID (4008,0210) is an example of a Data [Element](#) having Defined Terms. It is defined to have a [Value](#) that may be one of the set of standard Values; REPORT or AMENDMENT (see PS 3.3). Because this Data [Element](#) has Defined Terms other Interpretation [Type](#) IDs may be defined by the implementor.

```
#include <gdcmDefinedTerms.h>
```

### Public Member Functions

- [DefinedTerms](#) ()

### 27.77.1 Detailed Description

Defined Terms are used when the specified explicit Values may be extended by implementors to include additional new Values. These new Values shall be specified in the Conformance Statement (see PS 3.2) and shall not have the same meaning as currently defined Values in this standard. A Data [Element](#) with Defined Terms that does not contain a [Value](#) equivalent to one of the Values currently specified in this standard shall not be considered to have an invalid value. Note: Interpretation [Type](#) ID (4008,0210) is an example of a Data [Element](#) having Defined Terms. It is defined to have a [Value](#) that may be one of the set of standard Values; REPORT or AMENDMENT (see PS 3.3). Because this Data [Element](#) has Defined Terms other Interpretation [Type](#) IDs may be defined by the implementor.

### 27.77.2 Constructor & Destructor Documentation

27.77.2.1 `gdcm::DefinedTerms::DefinedTerms ( )` `[inline]`

The documentation for this class was generated from the following file:

- [gdcmDefinedTerms.h](#)

## 27.78 gdcm::Defs Class Reference

FIXME I do not like the name 'Defs'.

```
#include <gdcmDefs.h>
```

### Public Member Functions

- [Defs](#) ()
- [~Defs](#) ()
- const [IOD](#) & [GetIODFromFile](#) (const [File](#) &file) const
- const [IODs](#) & [GetIODs](#) () const
- [IODs](#) & [GetIODs](#) ()
- const [Macros](#) & [GetMacros](#) () const
- [Macros](#) & [GetMacros](#) ()

- const [Modules](#) & [GetModules](#) () const
- [Modules](#) & [GetModules](#) ()
- [Type](#) [GetTypeFromTag](#) (const [File](#) &file, const [Tag](#) &tag) const
- bool [IsEmpty](#) () const
- bool [Verify](#) (const [File](#) &file) const
- bool [Verify](#) (const [DataSet](#) &ds) const

### Static Public Member Functions

- static const char \* [GetIODNameFromMediaStorage](#) ([MediaStorage](#) const &ms)

### Protected Member Functions

- void [LoadDefaults](#) ()
- void [LoadFromFile](#) (const char \*filename)

### Friends

- class [Global](#)

## 27.78.1 Detailed Description

FIXME I do not like the name '[Defs](#)'.

#### Note

bla

#### Examples:

[GenerateStandardSOPClasses.cxx](#), and [TraverseModules.cxx](#).

## 27.78.2 Constructor & Destructor Documentation

27.78.2.1 `gdcm::Defs::Defs ( )`

27.78.2.2 `gdcm::Defs::~~Defs ( )`

## 27.78.3 Member Function Documentation

27.78.3.1 `const IOD& gdcm::Defs::GetIODFromFile ( const File & file ) const`

27.78.3.2 `static const char* gdcm::Defs::GetIODNameFromMediaStorage ( MediaStorage const & ms ) [static]`

#### Examples:

[GenerateStandardSOPClasses.cxx](#).

27.78.3.3 `const IODs& gdcM::Defs::GetIODs ( ) const [inline]`

Examples:

[TraverseModules.cxx](#).

27.78.3.4 `IODs& gdcM::Defs::GetIODs ( ) [inline]`

27.78.3.5 `const Macros& gdcM::Defs::GetMacros ( ) const [inline]`

Users should not directly use [Macro](#). [Macro](#) are simply a way for DICOM WG to re-use Tables. [Macros](#) are conveniently wrapped within [Modules](#). See [gdcM::Module](#) API directly

Examples:

[TraverseModules.cxx](#).

27.78.3.6 `Macros& gdcM::Defs::GetMacros ( ) [inline]`

27.78.3.7 `const Modules& gdcM::Defs::GetModules ( ) const [inline]`

Examples:

[TraverseModules.cxx](#).

27.78.3.8 `Modules& gdcM::Defs::GetModules ( ) [inline]`

27.78.3.9 `Type gdcM::Defs::GetTypeFromTag ( const File & file, const Tag & tag ) const`

27.78.3.10 `bool gdcM::Defs::IsEmpty ( ) const [inline]`

27.78.3.11 `void gdcM::Defs::LoadDefaults ( ) [protected]`

27.78.3.12 `void gdcM::Defs::LoadFromFile ( const char * filename ) [protected]`

27.78.3.13 `bool gdcM::Defs::Verify ( const File & file ) const`

27.78.3.14 `bool gdcM::Defs::Verify ( const DataSet & ds ) const`

## 27.78.4 Friends And Related Function Documentation

27.78.4.1 `friend class Global [friend]`

The documentation for this class was generated from the following file:

- [gdcMDefs.h](#)




## 27.79 gdcm::DeltaEncodingCodec Class Reference


[DeltaEncodingCodec](#) compression used by some private vendor.

```
#include <gdcmDeltaEncodingCodec.h>
```

Inheritance diagram for gdcm::DeltaEncodingCodec:

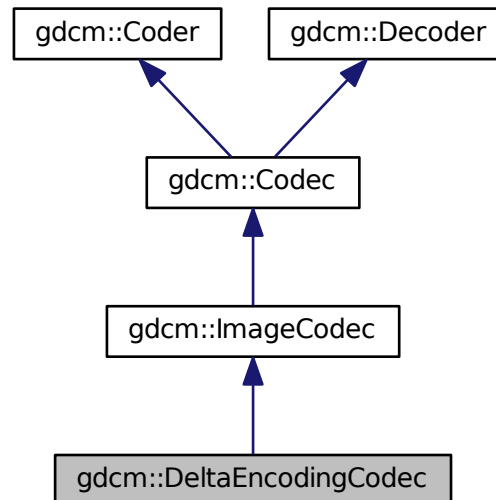

Collaboration diagram for gdcm::DeltaEncodingCodec:

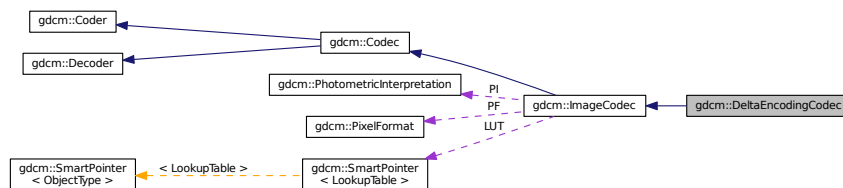

### Public Member Functions

- [DeltaEncodingCodec](#) ()
- [~DeltaEncodingCodec](#) ()
- bool [CanDecode](#) ([TransferSyntax](#) const &ts)
- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)

*Decode.*

## Protected Member Functions

- bool [Decode](#) (std::istream &is, std::ostream &os)

## Additional Inherited Members

### 27.79.1 Detailed Description

[DeltaEncodingCodec](#) compression used by some private vendor.

### 27.79.2 Constructor & Destructor Documentation

27.79.2.1 `gdcm::DeltaEncodingCodec::DeltaEncodingCodec ( )`

27.79.2.2 `gdcm::DeltaEncodingCodec::~~DeltaEncodingCodec ( )`

### 27.79.3 Member Function Documentation

27.79.3.1 `bool gdcm::DeltaEncodingCodec::CanDecode ( TransferSyntax const & ts )`

27.79.3.2 `bool gdcm::DeltaEncodingCodec::Decode ( DataElement const & , DataElement & )` [virtual]

Decode.

Reimplemented from [gdcm::Decoder](#).

27.79.3.3 `bool gdcm::DeltaEncodingCodec::Decode ( std::istream & is, std::ostream & os )` [protected]

The documentation for this class was generated from the following file:

- [gdcmDeltaEncodingCodec.h](#)

## 27.80 gdcm::DICOMDIR Class Reference

[DICOMDIR](#) class.

```
#include <gdcmDICOMDIR.h>
```

## Public Member Functions

- [DICOMDIR](#) ()
- [DICOMDIR](#) (const [FileSet](#) &fs)

### 27.80.1 Detailed Description

[DICOMDIR](#) class.

Structured for handling [DICOMDIR](#)

## 27.80.2 Constructor & Destructor Documentation

27.80.2.1 `gdcm::DICOMDIR::DICOMDIR ( )` [\[inline\]](#)

27.80.2.2 `gdcm::DICOMDIR::DICOMDIR ( const FileSet & fs )` [\[inline\]](#)

The documentation for this class was generated from the following file:

- [gdcmDICOMDIR.h](#)

## 27.81 gdcm::DICOMDIRGenerator Class Reference

[DICOMDIRGenerator](#) class This is a STD-GEN-CD [DICOMDIR](#) generator. ref: PS 3.11-2008 Annex D (Normative) - General Purpose CD-R and DVD Interchange Profiles.

```
#include <gdcmDICOMDIRGenerator.h>
```

### Public Types

- typedef [Directory::FileNamesType](#) FileNamesType
- typedef [Directory::FilenameType](#) FilenameType

### Public Member Functions

- [DICOMDIRGenerator](#) ()
- [~DICOMDIRGenerator](#) ()
- bool [Generate](#) ()  
*Main function to generate the [DICOMDIR](#).*
- [File](#) & [GetFile](#) ()
- void [SetDescriptor](#) (const char \*d)
- void [SetFile](#) (const [File](#) &f)  
*Set/Get file. The [DICOMDIR](#) file will be valid once a call to Generate has been done.*
- void [SetFileNames](#) ([FileNamesType](#) const &fns)  
*Set the list of filenames from which the [DICOMDIR](#) should be generated from.*
- void [SetRootDirectory](#) ([FilenameType](#) const &root)  
*Set the root directory from which the filenames should be considered.*

### Protected Member Functions

- bool [AddImageDirectoryRecord](#) ()
- bool [AddPatientDirectoryRecord](#) ()
- bool [AddSeriesDirectoryRecord](#) ()
- bool [AddStudyDirectoryRecord](#) ()
- [Scanner](#) & [GetScanner](#) ()

### 27.81.1 Detailed Description

[DICOMDIRGenerator](#) class This is a STD-GEN-CD [DICOMDIR](#) generator. ref: PS 3.11-2008 Annex D (Normative) - General Purpose CD-R and DVD Interchange Profiles.

#### Note

PS 3.11 - 2008 / D.3.2 Physical Medium And Medium Format The STD-GEN-CD and STD-GEN-SEC-CD application profiles require the 120 mm CD-R physical medium with the ISO/IEC 9660 Media Format, as defined in PS3.12. See also PS 3.12 - 2008 / Annex F 120mm CD-R Medium (Normative) and PS 3.10 - 2008 / 8 DICOM [File Service](#) / 8.1 FILE-SET

#### Warning

: PS 3.11 - 2008 / D.3.1 SOP Classes and Transfer Syntaxes Composite [Image](#) & Stand-alone Storage are required to be stored as Explicit [VR](#) Little Endian Uncompressed (1.2.840.10008.1.2.1). When a DICOM file is found using another Transfer Syntax the generator will simply stops.

- Input files should be Explicit [VR](#) Little Endian
- filenames should be valid [VR::CS](#) value (16 bytes, upper case ...)

**Bug** : There is a current limitation of not handling Referenced SOP Class UID / Referenced SOP Instance UID simply because the [Scanner](#) does not allow us See PS 3.11 / [Table D.3-2](#) STD-GEN Additional [DICOMDIR](#) Keys

### 27.81.2 Member Typedef Documentation

27.81.2.1 typedef [Directory::FilenameType](#) gdcm::DICOMDIRGenerator::FilenameType

27.81.2.2 typedef [Directory::FilenameType](#) gdcm::DICOMDIRGenerator::FilenameType

### 27.81.3 Constructor & Destructor Documentation

27.81.3.1 gdcm::DICOMDIRGenerator::DICOMDIRGenerator ( )

27.81.3.2 gdcm::DICOMDIRGenerator::~~DICOMDIRGenerator ( )

### 27.81.4 Member Function Documentation

27.81.4.1 bool gdcm::DICOMDIRGenerator::AddImageDirectoryRecord ( ) [protected]

27.81.4.2 bool gdcm::DICOMDIRGenerator::AddPatientDirectoryRecord ( ) [protected]

27.81.4.3 bool gdcm::DICOMDIRGenerator::AddSeriesDirectoryRecord ( ) [protected]

27.81.4.4 bool gdcm::DICOMDIRGenerator::AddStudyDirectoryRecord ( ) [protected]

27.81.4.5 bool gdcm::DICOMDIRGenerator::Generate ( )

Main function to generate the [DICOMDIR](#).

27.81.4.6 **File** & gdcm::DICOMDIRGenerator::GetFile ( )

27.81.4.7 **Scanner** & gdcm::DICOMDIRGenerator::GetScanner ( ) [protected]

27.81.4.8 void gdcm::DICOMDIRGenerator::SetDescriptor ( const char \* *d* )

Set the [File](#) Set ID.

Warning

this need to be a valid [VR::CS](#) value

27.81.4.9 void gdcm::DICOMDIRGenerator::SetFile ( const File & *f* )

Set/Get file. The [DICOMDIR](#) file will be valid once a call to Generate has been done.

27.81.4.10 void gdcm::DICOMDIRGenerator::SetFilenames ( [FilenameType](#) const & *fns* )

Set the list of filenames from which the [DICOMDIR](#) should be generated from.

27.81.4.11 void gdcm::DICOMDIRGenerator::SetRootDirectory ( [FilenameType](#) const & *root* )

Set the root directory from which the filenames should be considered.

The documentation for this class was generated from the following file:

- [gdcmDICOMDIRGenerator.h](#)

## 27.82 gdcm::Dict Class Reference

Class to represent a map of [DictEntry](#).

```
#include <gdcmDict.h>
```

### Public Types

- typedef MapDictEntry::const\_iterator [ConstIterator](#)
- typedef MapDictEntry::iterator [Iterator](#)
- typedef std::map< [Tag](#), [DictEntry](#) > [MapDictEntry](#)

### Public Member Functions

- [Dict](#) ()
- void [AddDictEntry](#) (const [Tag](#) &tag, const [DictEntry](#) &de)
- [ConstIterator](#) [Begin](#) () const
- [ConstIterator](#) [End](#) () const
- const [DictEntry](#) & [GetDictEntry](#) (const [Tag](#) &tag) const
- const [DictEntry](#) & [GetDictEntryByKeyword](#) (const char \*keyword, [Tag](#) &tag) const

- const [DictEntry](#) & [GetDictEntryByName](#) (const char \*name, [Tag](#) &tag) const
- const char \* [GetKeywordFromTag](#) ([Tag](#) const &tag) const  
*Function to return the Keyword from a [Tag](#).*
- bool [IsEmpty](#) () const

## Protected Member Functions

- void [LoadDefault](#) ()

## Friends

- class [Dicts](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [Dict](#) &\_val)

### 27.82.1 Detailed Description

Class to represent a map of [DictEntry](#).

#### Note

bla TODO FIXME: For [Element](#) == 0x0 need to return Name = Group Length ValueRepresentation = UL Value← Multiplicity = 1

#### Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [PublicDict.cxx](#), and [ReadAndPrintAttributes.cxx](#).

### 27.82.2 Member Typedef Documentation

27.82.2.1 typedef MapDictEntry::const\_iterator [gdcm::Dict::ConstIterator](#)

27.82.2.2 typedef MapDictEntry::iterator [gdcm::Dict::Iterator](#)

27.82.2.3 typedef std::map<[Tag](#), [DictEntry](#)> [gdcm::Dict::MapDictEntry](#)

### 27.82.3 Constructor & Destructor Documentation

27.82.3.1 [gdcm::Dict::Dict](#) ( ) [\[inline\]](#)

### 27.82.4 Member Function Documentation

27.82.4.1 void [gdcm::Dict::AddDictEntry](#) ( const [Tag](#) & tag, const [DictEntry](#) & de ) [\[inline\]](#)

27.82.4.2 [ConstIterator](#) [gdcm::Dict::Begin](#) ( ) const [\[inline\]](#)

#### Examples:

[GenAllVR.cxx](#), and [GenFakeIdentifyFile.cxx](#).

27.82.4.3 `ConstIterator gdcmm::Dict::End ( ) const` `[inline]`

Examples:

[GenAllVR.cxx](#), and [GenFakeIdentifyFile.cxx](#).

27.82.4.4 `const DictEntry& gdcmm::Dict::GetDictEntry ( const Tag & tag ) const` `[inline]`

Examples:

[GenFakeIdentifyFile.cxx](#), and [PublicDict.cxx](#).

27.82.4.5 `const DictEntry& gdcmm::Dict::GetDictEntryByKeyword ( const char * keyword, Tag & tag ) const` `[inline]`

Lookup [DictEntry](#) by keyword. Even if DICOM standard defines keyword as being unique. The lookup table is built on [Tag](#). Therefore looking up a [DictEntry](#) by Keyword is more inefficient than looking up by [Tag](#).

27.82.4.6 `const DictEntry& gdcmm::Dict::GetDictEntryByName ( const char * name, Tag & tag ) const` `[inline]`

Inefficient way of looking up tag by name. Technically DICOM does not guarantee uniqueness (and [Curve](#) / [Overlay](#) are there to prove it). But most of the time name is in fact unique and can be uniquely link to a tag

Examples:

[ReadAndPrintAttributes.cxx](#).

27.82.4.7 `const char* gdcmm::Dict::GetKeywordFromTag ( Tag const & tag ) const` `[inline]`

Function to return the Keyword from a [Tag](#).

27.82.4.8 `bool gdcmm::Dict::IsEmpty ( ) const` `[inline]`

Referenced by `gdcmm::Dicts::IsEmpty()`.

27.82.4.9 `void gdcmm::Dict::LoadDefault ( )` `[protected]`

## 27.82.5 Friends And Related Function Documentation

27.82.5.1 `friend class Dicts` `[friend]`

27.82.5.2 `std::ostream& operator<< ( std::ostream & _os, const Dict & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmmDict.h](#)

## 27.83 gdcmmDictConverter Class Reference

Class to convert a .dic file into something else:

```
#include <gdcmmDictConverter.h>
```

### Public Types

- enum [OutputTypes](#) {  
    [DICT\\_DEFAULT](#) = 0,  
    [DICT\\_DEBUG](#),  
    [DICT\\_XML](#) }

### Public Member Functions

- [DictConverter](#) ()
- [~DictConverter](#) ()
- void [Convert](#) ()
- const std::string & [GetDictName](#) () const
- const std::string & [GetInputFilename](#) () const
- const std::string & [GetOutputFilename](#) () const
- int [GetOutputType](#) () const
- void [SetDictName](#) (const char \*name)
- void [SetInputFileName](#) (const char \*filename)
- void [SetOutputFileName](#) (const char \*filename)
- void [SetOutputType](#) (int type)

### Static Public Member Functions

- static bool [Readuint16](#) (const char \*raw, uint16\_t &ov)
- static bool [ReadVM](#) (const char \*raw, [VM::VMType](#) &type)
- static bool [ReadVR](#) (const char \*raw, [VR::VRType](#) &type)

### Protected Member Functions

- void [AddGroupLength](#) ()
- bool [ConvertToCXX](#) (const char \*raw, std::string &cxx)
- bool [ConvertToXML](#) (const char \*raw, std::string &cxx)
- void [WriteFooter](#) ()
- void [WriteHeader](#) ()

#### 27.83.1 Detailed Description

Class to convert a .dic file into something else:

- CXX code : embeded dict into shared lib ([DICT\\_DEFAULT](#))
- Debug mode ([DICT\\_DEBUG](#))
- XML dict ([DICT\\_XML](#))



Note

## 27.83.2 Member Enumeration Documentation

### 27.83.2.1 enum gdcm::DictConverter::OutputTypes

Enumerator

***DICT\_DEFAULT***

***DICT\_DEBUG***

***DICT\_XML***

## 27.83.3 Constructor & Destructor Documentation

27.83.3.1 gdcm::DictConverter::DictConverter ( )

27.83.3.2 gdcm::DictConverter::~~DictConverter ( )

## 27.83.4 Member Function Documentation

27.83.4.1 void gdcm::DictConverter::AddGroupLength ( ) [protected]

27.83.4.2 void gdcm::DictConverter::Convert ( )

27.83.4.3 bool gdcm::DictConverter::ConvertToCXX ( const char \* raw, std::string & cxx ) [protected]

27.83.4.4 bool gdcm::DictConverter::ConvertToXML ( const char \* raw, std::string & cxx ) [protected]

27.83.4.5 const std::string& gdcm::DictConverter::GetDictName ( ) const

27.83.4.6 const std::string& gdcm::DictConverter::GetInputFilename ( ) const

27.83.4.7 const std::string& gdcm::DictConverter::GetOutputFilename ( ) const

27.83.4.8 int gdcm::DictConverter::GetOutputType ( ) const [inline]

27.83.4.9 static bool gdcm::DictConverter::Readuint16 ( const char \* raw, uint16\_t & ov ) [static]

27.83.4.10 static bool gdcm::DictConverter::ReadVM ( const char \* raw, VM::VMType & type ) [static]

27.83.4.11 static bool gdcm::DictConverter::ReadVR ( const char \* raw, VR::VRType & type ) [static]

27.83.4.12 void gdcm::DictConverter::SetDictName ( const char \* name )

27.83.4.13 void gdcm::DictConverter::SetInputFileName ( const char \* filename )

27.83.4.14 void gdcm::DictConverter::SetOutputFileName ( const char \* filename )

27.83.4.15 void gdcm::DictConverter::SetOutputType ( int type ) [inline]

27.83.4.16 void gdcmmDictConverter::WriteFooter ( ) [protected]

27.83.4.17 void gdcmmDictConverter::WriteHeader ( ) [protected]

The documentation for this class was generated from the following file:

- [gdcmmDictConverter.h](#)

## 27.84 gdcmmDictEntry Class Reference

Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcmmTag](#) to the needed information.

```
#include <gdcmmDictEntry.h>
```

### Public Member Functions

- [DictEntry](#) (const char \*name="", const char \*keyword="", [VR](#) const &vr=[VR::INVALID](#), [VM](#) const &vm=[VM::VM0](#), bool ret=false)
- const char \* [GetKeyword](#) () const  
*same as GetName but without spaces...*
- const char \* [GetName](#) () const  
*Set/Get Name.*
- bool [GetRetired](#) () const  
*Set/Get Retired flag.*
- const [VM](#) & [GetVM](#) () const  
*Set/Get VM.*
- const [VR](#) & [GetVR](#) () const  
*Set/Get VR.*
- bool [IsUnique](#) () const
- void [SetElementXX](#) (bool v)  
*Set whether element is shared in multiple elements (Source [Image](#) IDs typically)*
- void [SetGroupXX](#) (bool v)  
*Set whether element is shared in multiple groups (Curve/Overlay typically)*
- void [SetKeyword](#) (const char \*keyword)
- void [SetName](#) (const char \*name)
- void [SetRetired](#) (bool retired)
- void [SetVM](#) ([VM](#) const &vm)
- void [SetVR](#) (const [VR](#) &vr)

### Friends

- class [Dict](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [DictEntry](#) &\_val)

### 27.84.1 Detailed Description

Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcmm::Tag](#) to the needed information.

#### Note

bla TODO FIXME: Need a PublicDictEntry...indeed [DictEntry](#) has a notion of retired which does not exist in PrivateDictEntry...

#### See also

[gdcmm::Dict](#)

#### Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [PublicDict.cxx](#), and [TraverseModules.cxx](#).

### 27.84.2 Constructor & Destructor Documentation

27.84.2.1 `gdcmm::DictEntry::DictEntry ( const char * name = " ", const char * keyword = " ", VR const & vr = VR::INVALID, VM const & vm = VM::VMO, bool ret = false ) [inline]`

### 27.84.3 Member Function Documentation

27.84.3.1 `const char* gdcmm::DictEntry::GetKeyword ( ) const [inline]`

same as GetName but without spaces...

27.84.3.2 `const char* gdcmm::DictEntry::GetName ( ) const [inline]`

Set/Get Name.

Referenced by `gdcmm::PrivateDict::PrintXML()`.

27.84.3.3 `bool gdcmm::DictEntry::GetRetired ( ) const [inline]`

Set/Get Retired flag.

#### Examples:

[GenAllVR.cxx](#).

27.84.3.4 `const VM& gdcmm::DictEntry::GetVM ( ) const [inline]`

Set/Get [VM](#).

Referenced by `gdcmm::PrivateDict::AddDictEntry()`, and `gdcmm::PrivateDict::PrintXML()`.

27.84.3.5 `const VR& gdcmm::DictEntry::GetVR ( ) const` `[inline]`

Set/Get [VR](#).

Examples:

[GenAllVR.cxx](#), and [GenFakeIdentifyFile.cxx](#).

Referenced by `gdcmm::PrivateDict::AddDictEntry()`, and `gdcmm::PrivateDict::PrintXML()`.

27.84.3.6 `bool gdcmm::DictEntry::IsUnique ( ) const` `[inline]`

Return whether the name of the [DataElement](#) can be considered to be unique. As of 2008 all elements name were unique (except the explicitly 'XX' ones)

27.84.3.7 `void gdcmm::DictEntry::SetElementXX ( bool v )` `[inline]`

Set whether element is shared in multiple elements (Source [Image](#) IDs typically)

27.84.3.8 `void gdcmm::DictEntry::SetGroupXX ( bool v )` `[inline]`

Set whether element is shared in multiple groups (Curve/Overlay typically)

27.84.3.9 `void gdcmm::DictEntry::SetKeyword ( const char * keyword )` `[inline]`

27.84.3.10 `void gdcmm::DictEntry::SetName ( const char * name )` `[inline]`

27.84.3.11 `void gdcmm::DictEntry::SetRetired ( bool retired )` `[inline]`

27.84.3.12 `void gdcmm::DictEntry::SetVM ( VM const & vm )` `[inline]`

27.84.3.13 `void gdcmm::DictEntry::SetVR ( const VR & vr )` `[inline]`

Referenced by `gdcmm::PrivateDict::AddDictEntry()`.

## 27.84.4 Friends And Related Function Documentation

27.84.4.1 `friend class Dict` `[friend]`

27.84.4.2 `std::ostream& operator<< ( std::ostream &_os, const DictEntry &_val )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmmDictEntry.h](#)

## 27.85 gdcmm::DictPrinter Class Reference


[DictPrinter](#) class.

```
#include <gdcmmDictPrinter.h>
```

Inheritance diagram for gdcmmDictPrinter:

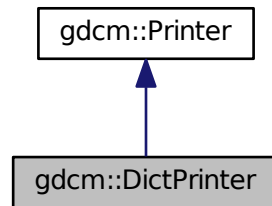

Collaboration diagram for gdcmmDictPrinter:

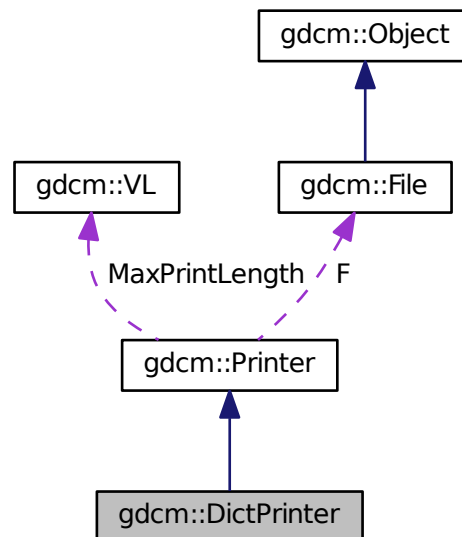

## Public Member Functions

- [DictPrinter](#) ()
- [~DictPrinter](#) ()
- void [Print](#) (std::ostream &os)

## Protected Member Functions

- void [PrintDataElement2](#) (std::ostream &os, const [DataSet](#) &ds, const [DataElement](#) &ide)
- void [PrintDataSet2](#) (std::ostream &os, const [DataSet](#) &ds)

## Additional Inherited Members

### 27.85.1 Detailed Description

[DictPrinter](#) class.

### 27.85.2 Constructor & Destructor Documentation

27.85.2.1 `gdcmm::DictPrinter::DictPrinter ( )`

27.85.2.2 `gdcmm::DictPrinter::~~DictPrinter ( )`

### 27.85.3 Member Function Documentation

27.85.3.1 `void gdcmm::DictPrinter::Print ( std::ostream & os )`

27.85.3.2 `void gdcmm::DictPrinter::PrintDataElement2 ( std::ostream & os, const DataSet & ds, const DataElement & ide )`  
[protected]

27.85.3.3 `void gdcmm::DictPrinter::PrintDataSet2 ( std::ostream & os, const DataSet & ds )` [protected]

The documentation for this class was generated from the following file:

- [gdcmmDictPrinter.h](#)

## 27.86 gdcmm::Dicts Class Reference

Class to manipulate the sum of knowledge (all the dict user load)

```
#include <gdcmmDicts.h>
```

## Public Member Functions

- [Dicts](#) ()
- [~Dicts](#) ()
- const [CSAHeaderDict](#) & [GetCSAHeaderDict](#) () const
- const [DictEntry](#) & [GetDictEntry](#) (const [Tag](#) &tag, const char \*owner=NULL) const  
*NOT THREAD SAFE.*
- const [DictEntry](#) & [GetDictEntry](#) (const [PrivateTag](#) &tag) const
- const [PrivateDict](#) & [GetPrivateDict](#) () const
- [PrivateDict](#) & [GetPrivateDict](#) ()
- const [Dict](#) & [GetPublicDict](#) () const
- bool [IsEmpty](#) () const

## Protected Types

- enum [ConstructorType](#) {  
    [PHILIPS](#),  
    [GEMS](#),  
    [SIEMENS](#) }

## Protected Member Functions

- void [LoadDefaults](#) ()

## Static Protected Member Functions

- static const char \* [GetConstructorString](#) ([ConstructorType](#) type)

## Friends

- class [Global](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [Dicts](#) &d)

### 27.86.1 Detailed Description

Class to manipulate the sum of knowledge (all the dict user load)

#### Note

bla

#### Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [PublicDict.cxx](#), [ReadAndPrintAttributes.cxx](#), and [TraverseModules.cxx](#).

### 27.86.2 Member Enumeration Documentation

27.86.2.1 enum `gdcm::Dicts::ConstructorType` [protected]

#### Enumerator

***PHILIPS***

***GEMS***

***SIEMENS***

### 27.86.3 Constructor & Destructor Documentation

27.86.3.1 `gdcm::Dicts::Dicts ( )`

27.86.3.2 `gdcm::Dicts::~~Dicts ( )`

### 27.86.4 Member Function Documentation

27.86.4.1 `static const char* gdcmm::Dicts::GetConstructorString ( ConstructorType type )` `[static],[protected]`

27.86.4.2 `const CSAHeaderDict& gdcmm::Dicts::GetCSAHeaderDict ( ) const`

Examples:

[MrProtocol.cxx](#).

27.86.4.3 `const DictEntry& gdcmm::Dicts::GetDictEntry ( const Tag & tag, const char * owner = NULL ) const`

NOT THREAD SAFE.

works for both public and private dicts: owner is null for public dict

Warning

owner need to be set to appropriate owner for call to work. see

Examples:

[PublicDict.cxx](#), and [TraverseModules.cxx](#).

27.86.4.4 `const DictEntry& gdcmm::Dicts::GetDictEntry ( const PrivateTag & tag ) const`

27.86.4.5 `const PrivateDict& gdcmm::Dicts::GetPrivateDict ( ) const`

27.86.4.6 `PrivateDict& gdcmm::Dicts::GetPrivateDict ( )`

27.86.4.7 `const Dict& gdcmm::Dicts::GetPublicDict ( ) const`

Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [PublicDict.cxx](#), and [ReadAndPrintAttributes.cxx](#).

27.86.4.8 `bool gdcmm::Dicts::IsEmpty ( ) const` `[inline]`

References `gdcmm::Dict::IsEmpty()`.

27.86.4.9 `void gdcmm::Dicts::LoadDefaults ( )` `[protected]`

## 27.86.5 Friends And Related Function Documentation

27.86.5.1 `friend class Global` `[friend]`

27.86.5.2 `std::ostream& operator<< ( std::ostream & _os, const Dicts & d )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmmDicts.h](#)



## 27.87 gdcm::network::DIMSE Class Reference

DIMSE PS 3.7 - 2009 Annex E [Command](#) Dictionary (Normative) E.1 REGISTRY OF DICOM COMMAND ELEMENTS  
[Table E.1-1](#) COMMAND FIELDS (PART 1)

```
#include <gdcmDIMSE.h>
```

### Public Types

- enum [CommandTypes](#) {  
    [C\\_STORE\\_RQ](#) = 0x0001,  
    [C\\_STORE\\_RSP](#) = 0x8001,  
    [C\\_GET\\_RQ](#) = 0x0010,  
    [C\\_GET\\_RSP](#) = 0x8010,  
    [C\\_FIND\\_RQ](#) = 0x0020,  
    [C\\_FIND\\_RSP](#) = 0x8020,  
    [C\\_MOVE\\_RQ](#) = 0x0021,  
    [C\\_MOVE\\_RSP](#) = 0x8021,  
    [C\\_ECHO\\_RQ](#) = 0x0030,  
    [C\\_ECHO\\_RSP](#) = 0x8030,  
    [N\\_EVENT\\_REPORT\\_RQ](#) = 0x0100,  
    [N\\_EVENT\\_REPORT\\_RSP](#) = 0x8100,  
    [N\\_GET\\_RQ](#) = 0x0110,  
    [N\\_GET\\_RSP](#) = 0x8110,  
    [N\\_SET\\_RQ](#) = 0x0120,  
    [N\\_SET\\_RSP](#) = 0x8120,  
    [N\\_ACTION\\_RQ](#) = 0x0130,  
    [N\\_ACTION\\_RSP](#) = 0x8130,  
    [N\\_CREATE\\_RQ](#) = 0x0140,  
    [N\\_CREATE\\_RSP](#) = 0x8140,  
    [N\\_DELETE\\_RQ](#) = 0x0150,  
    [N\\_DELETE\\_RSP](#) = 0x8150,  
    [C\\_CANCEL\\_RQ](#) = 0x0FFF }

### 27.87.1 Detailed Description

DIMSE PS 3.7 - 2009 Annex E [Command](#) Dictionary (Normative) E.1 REGISTRY OF DICOM COMMAND ELEMENTS  
[Table E.1-1](#) COMMAND FIELDS (PART 1)

### 27.87.2 Member Enumeration Documentation

#### 27.87.2.1 enum gdcm::network::DIMSE::CommandTypes

Enumerator

***C\_STORE\_RQ***  
***C\_STORE\_RSP***  
***C\_GET\_RQ***  
***C\_GET\_RSP***  
***C\_FIND\_RQ***  
***C\_FIND\_RSP***

**C\_MOVE\_RQ**  
**C\_MOVE\_RSP**  
**C\_ECHO\_RQ**  
**C\_ECHO\_RSP**  
**N\_EVENT\_REPORT\_RQ**  
**N\_EVENT\_REPORT\_RSP**  
**N\_GET\_RQ**  
**N\_GET\_RSP**  
**N\_SET\_RQ**  
**N\_SET\_RSP**  
**N\_ACTION\_RQ**  
**N\_ACTION\_RSP**  
**N\_CREATE\_RQ**  
**N\_CREATE\_RSP**  
**N\_DELETE\_RQ**  
**N\_DELETE\_RSP**  
**C\_CANCEL\_RQ**

The documentation for this class was generated from the following file:

- [gdcmdIMSE.h](#)

## 27.88 gdcmd::DirectionCosines Class Reference

class to handle [DirectionCosines](#)

```
#include <gdcmdDirectionCosines.h>
```

### Public Member Functions

- [DirectionCosines](#) ()
- [DirectionCosines](#) (const double dircos[6])
- [~DirectionCosines](#) ()
- double [ComputeDistAlongNormal](#) (const double ipp[3]) const  
*Compute the distance along the normal.*
- void [Cross](#) (double z[3]) const  
*Compute Cross product.*
- double [CrossDot](#) ([DirectionCosines](#) const &dc) const  
*Compute the Dot product of the two cross vector of both [DirectionCosines](#) object.*
- double [Dot](#) () const  
*Compute Dot.*
- bool [IsValid](#) () const  
*Return whether or not this is a valid direction cosines.*
- void [Normalize](#) ()  
*Normalize in-place.*

- `operator const double * () const`  
*Make the class behave like a const double \*.*
- `void Print (std::ostream &) const`  
*Print.*
- `bool SetFromString (const char *str)`

### 27.88.1 Detailed Description

class to handle [DirectionCosines](#)

Examples:

[DiscriminateVolume.cxx](#).

### 27.88.2 Constructor & Destructor Documentation

27.88.2.1 `gdcm::DirectionCosines::DirectionCosines ( )`

27.88.2.2 `gdcm::DirectionCosines::DirectionCosines ( const double dircos[6] )`

27.88.2.3 `gdcm::DirectionCosines::~~DirectionCosines ( )`

### 27.88.3 Member Function Documentation

27.88.3.1 `double gdcm::DirectionCosines::ComputeDistAlongNormal ( const double ipp[3] ) const`

Compute the distance along the normal.

27.88.3.2 `void gdcm::DirectionCosines::Cross ( double z[3] ) const`

Compute Cross product.

27.88.3.3 `double gdcm::DirectionCosines::CrossDot ( DirectionCosines const & dc ) const`

Compute the Dot product of the two cross vector of both [DirectionCosines](#) object.

Examples:

[DiscriminateVolume.cxx](#).

27.88.3.4 `double gdcm::DirectionCosines::Dot ( ) const`

Compute Dot.

27.88.3.5 `bool gdcm::DirectionCosines::IsValid ( ) const`

Return whether or not this is a valid direction cosines.

27.88.3.6 void `gdcmm::DirectionCosines::Normalize ( )`

Normalize in-place.

27.88.3.7 `gdcmm::DirectionCosines::operator const double * ( ) const` `[inline]`

Make the class behave like a const double \*.

27.88.3.8 void `gdcmm::DirectionCosines::Print ( std::ostream & ) const`

Print.

27.88.3.9 bool `gdcmm::DirectionCosines::SetFromString ( const char * str )`

Initialize from string str. It requires 6 floating point separated by a backslash character.

Examples:

[DiscriminateVolume.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmmDirectionCosines.h](#)

## 27.89 gdcmm::Directory Class Reference

Class for manipulation directories.

```
#include <gdcmmDirectory.h>
```

### Public Types

- typedef std::vector< [FilenameType](#) > [FileNamesType](#)
- typedef std::string [FilenameType](#)

### Public Member Functions

- [Directory](#) ()
- [~Directory](#) ()
- [FileNamesType](#) const & [GetDirectories](#) () const  
*Return the Directories traversed.*
- [FileNamesType](#) const & [GetFileNames](#) () const  
*Set/Get the file names within the directory.*
- [FilenameType](#) const & [GetToplevel](#) () const  
*Get the name of the toplevel directory.*
- unsigned int [Load](#) ([FilenameType](#) const &name, bool recursive=false)
- void [Print](#) (std::ostream &os=std::cout) const  
*Print.*

## Protected Member Functions

- unsigned int [Explore](#) ([FilenameType](#) const &name, bool recursive)  
*Return number of file found when 'recursive'ly exploring directory `name`*

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Directory](#) &d)

### 27.89.1 Detailed Description

Class for manipulation directories.

#### Note

This implementation provide a cross platform implementation for manipulating diretores: basically traversing directories and harvesting files  
will not take into account unix type hidden file recursive option will not look into UNIX type hidden directory (those starting with a '.')  
Since python or C# provide there own equivalent implementation, in which case [gdcm::Directory](#) does not make much sense.

#### Examples:

[DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [gdcmorthoplanes.cxx](#), [GenerateRTSTRUCT.cxx](#), [ReadUTF8Qt←Dir.cxx](#), [reslicesphere.cxx](#), [SortImage.cxx](#), [threadgdcm.cxx](#), and [VolumeSorter.cxx](#).

### 27.89.2 Member Typedef Documentation

27.89.2.1 `typedef std::vector<FilenameType> gdcm::Directory::FilenamesType`

#### Examples:

[DiscriminateVolume.cxx](#).

27.89.2.2 `typedef std::string gdcm::Directory::FilenameType`

### 27.89.3 Constructor & Destructor Documentation

27.89.3.1 `gdcm::Directory::Directory ( )` `[inline]`

27.89.3.2 `gdcm::Directory::~~Directory ( )` `[inline]`

### 27.89.4 Member Function Documentation

27.89.4.1 `unsigned int gdcm::Directory::Explore ( FilenameType const & name, bool recursive )` `[protected]`

Return number of file found when 'recursive'ly exploring directory `name`

**27.89.4.2** `FilenameType const& gdcmm::Directory::GetDirectories ( ) const` `[inline]`

Return the Directories traversed.

**27.89.4.3** `FilenameType const& gdcmm::Directory::GetFileNames ( ) const` `[inline]`

Set/Get the file names within the directory.

Examples:

[DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [gdcmmorphoplanes.cxx](#), [GenerateRTSTRUCT.cxx](#), [ReadUTF8QtDir.cxx](#), [reslicesphere.cxx](#), [SortImage.cxx](#), [threadgdcmm.cxx](#), and [VolumeSorter.cxx](#).

**27.89.4.4** `FilenameType const& gdcmm::Directory::GetToplevel ( ) const` `[inline]`

Get the name of the toplevel directory.

**27.89.4.5** `unsigned int gdcmm::Directory::Load ( FilenameType const & name, bool recursive = false )` `[inline]`

construct a list of filenames and subdirectory beneath directory: name

Warning

: hidden file and hidden directory are not loaded.

Examples:

[DecompressImageMultiframe.cs](#), [DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [gdcmmorphoplanes.cxx](#), [GenerateRTSTRUCT.cxx](#), [ReadUTF8QtDir.cxx](#), [reslicesphere.cxx](#), [SortImage.cxx](#), [threadgdcmm.cxx](#), and [VolumeSorter.cxx](#).

**27.89.4.6** `void gdcmm::Directory::Print ( std::ostream & os = std::cout ) const`

Print.

Examples:

[SortImage.cxx](#).

Referenced by `gdcmm::operator<<()`.

## 27.89.5 Friends And Related Function Documentation

**27.89.5.1** `std::ostream& operator<< ( std::ostream & _os, const Directory & d )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmmDirectory.h](#)

## 27.90 gdcm::DirectoryHelper Class Reference

[DirectoryHelper](#) this class is designed to help mitigate some of the commonly performed operations on directories. namely: 1) the ability to determine the number of series in a directory by what type of series is present 2) the ability to find all ct series in a directory 3) the ability to find all mr series in a directory 4) to load a set of DataSets from a series that's already been sorted by the IPP sorter 5) For rtstruct stuff, you need to know the sopinstanceuid of each z plane, so there's a retrieval function for that 6) then a few other functions for rtstruct writeouts.

```
#include <gdcmDirectoryHelper.h>
```

### Static Public Member Functions

- static [Directory::FilenamesType GetCTImageSeriesUIDs](#) (const std::string &inDirectory)
- static [Directory::FilenamesType GetFilenamesFromSeriesUIDs](#) (const std::string &inDirectory, const std::string &inSeriesUID)
- static std::string [GetFrameOfReference](#) (const std::vector< [DataSet](#) > &inDS)
- static [Directory::FilenamesType GetMRIImageSeriesUIDs](#) (const std::string &inDirectory)
- static [Directory::FilenamesType GetRTStructSeriesUIDs](#) (const std::string &inDirectory)
- static [Directory::FilenamesType GetSeriesUIDsBySOPClassUID](#) (const std::string &inDirectory, const std::string &inSOPClassUID)
- static std::string [GetSOPClassUID](#) (const std::vector< [DataSet](#) > &inDS)
- static std::string [GetStringValueFromTag](#) (const [Tag](#) &t, const [DataSet](#) &ds)
- static std::vector< [DataSet](#) > [LoadImageFromFiles](#) (const std::string &inDirectory, const std::string &inSeriesUID)
- static std::string [RetrieveSOPInstanceUIDFromIndex](#) (int inIndex, const std::vector< [DataSet](#) > &inDS)
- static std::string [RetrieveSOPInstanceUIDFromZPosition](#) (double inZPos, const std::vector< [DataSet](#) > &inDS)

### 27.90.1 Detailed Description

[DirectoryHelper](#) this class is designed to help mitigate some of the commonly performed operations on directories. namely: 1) the ability to determine the number of series in a directory by what type of series is present 2) the ability to find all ct series in a directory 3) the ability to find all mr series in a directory 4) to load a set of DataSets from a series that's already been sorted by the IPP sorter 5) For rtstruct stuff, you need to know the sopinstanceuid of each z plane, so there's a retrieval function for that 6) then a few other functions for rtstruct writeouts.

### 27.90.2 Member Function Documentation

**27.90.2.1** static [Directory::FilenamesType](#) [gdcm::DirectoryHelper::GetCTImageSeriesUIDs](#) ( const std::string & *inDirectory* )  
[static]

**27.90.2.2** static [Directory::FilenamesType](#) [gdcm::DirectoryHelper::GetFilenamesFromSeriesUIDs](#) ( const std::string & *inDirectory*, const std::string & *inSeriesUID* ) [static]

Examples:

[GenerateRTSTRUCT.cxx](#).

**27.90.2.3** static std::string [gdcm::DirectoryHelper::GetFrameOfReference](#) ( const std::vector< [DataSet](#) > & *inDS* ) [static]

**27.90.2.4** static [Directory::FilenamesType](#) [gdcm::DirectoryHelper::GetMRIImageSeriesUIDs](#) ( const std::string & *inDirectory* )  
[static]

27.90.2.5 **static** **Directory::FilenameType** **gdcm::DirectoryHelper::GetRTStructSeriesUIDs** ( **const** **std::string** & *inDirectory* )  
[static]

Examples:

[GenerateRTSTRUCT.cxx](#).

27.90.2.6 **static** **Directory::FilenameType** **gdcm::DirectoryHelper::GetSeriesUIDsBySOPClassUID** ( **const** **std::string** & *inDirectory*, **const** **std::string** & *inSOPClassUID* ) [static]

27.90.2.7 **static** **std::string** **gdcm::DirectoryHelper::GetSOPClassUID** ( **const** **std::vector**< **DataSet** > & *inDS* ) [static]

27.90.2.8 **static** **std::string** **gdcm::DirectoryHelper::GetStringValueFromTag** ( **const** **Tag** & *t*, **const** **DataSet** & *ds* ) [static]

27.90.2.9 **static** **std::vector**<**DataSet**> **gdcm::DirectoryHelper::LoadImageFromFiles** ( **const** **std::string** & *inDirectory*, **const** **std::string** & *inSeriesUID* ) [static]

27.90.2.10 **static** **std::string** **gdcm::DirectoryHelper::RetrieveSOPInstanceUIDFromIndex** ( **int** *inIndex*, **const** **std::vector**< **DataSet** > & *inDS* ) [static]

27.90.2.11 **static** **std::string** **gdcm::DirectoryHelper::RetrieveSOPInstanceUIDFromZPosition** ( **double** *inZPos*, **const** **std::vector**< **DataSet** > & *inDS* ) [static]

The documentation for this class was generated from the following file:

- [gdcmDirectoryHelper.h](#)

## 27.91 gdcm::DummyValueGenerator Class Reference

Class for generating dummy value.

```
#include <gdcmDummyValueGenerator.h>
```

### Static Public Member Functions

- **static** **const** **char** \* [Generate](#) (**const** **char** \**input*)

#### 27.91.1 Detailed Description

Class for generating dummy value.

See also

[Anonymizer](#)

#### 27.91.2 Member Function Documentation



27.91.2.1 `static const char* gdcm::DummyValueGenerator::Generate ( const char * input )` `[static]`

Generate a dummy value from an input value. This is guarantee to always return the same output value when input is identical. Return an array of bytes that can be used for anonymization purpose, return NULL on error NOT THREAD SAFE

The documentation for this class was generated from the following file:


- [gdcmDummyValueGenerator.h](#)

## 27.92 gdcm::Dumper Class Reference


[Codec](#) class.

```
#include <gdcmDumper.h>
```

Inheritance diagram for `gdcm::Dumper`:

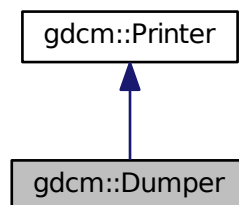

Collaboration diagram for `gdcM::Dumper`:

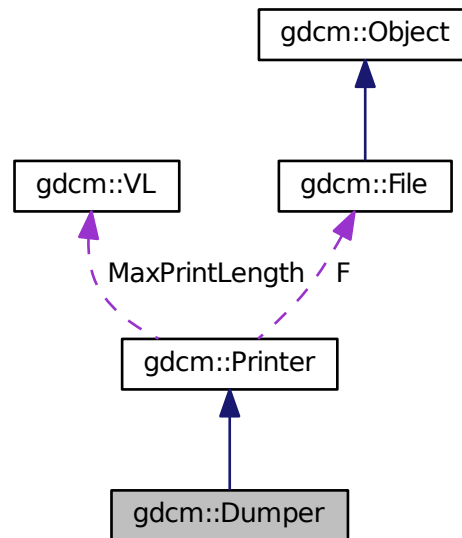

## Public Member Functions

- [Dumper \(\)](#)
- [~Dumper \(\)](#)

## Additional Inherited Members

### 27.92.1 Detailed Description

[Codec](#) class.

#### Note

Use it to simply dump value read from the file. No interpretation is done. But it is real fast ! Almost no overhead

### 27.92.2 Constructor & Destructor Documentation

27.92.2.1 `gdcM::Dumper::Dumper ( )` `[inline]`

27.92.2.2 `gdcM::Dumper::~~Dumper ( )` `[inline]`

The documentation for this class was generated from the following file:


- [gdcMDumper.h](#)

## 27.93 gdcmm::Element< TVR, TVM > Class Template Reference


[Element](#) class.

```
#include <gdcmmElement.h>
```

Inheritance diagram for gdcmm::Element< TVR, TVM >:

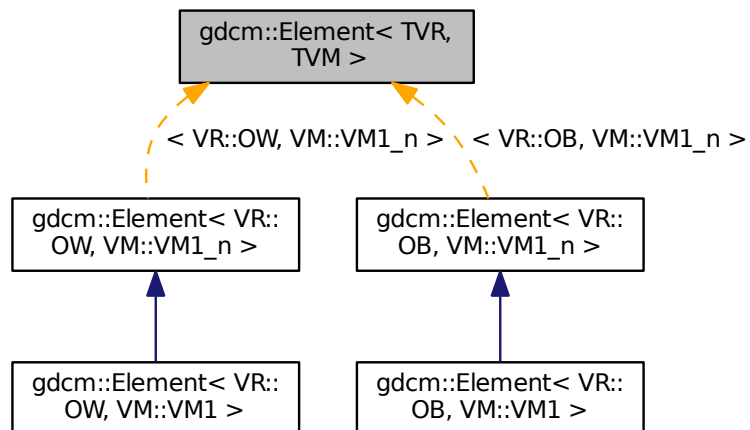

Collaboration diagram for gdcmm::Element< TVR, TVM >:

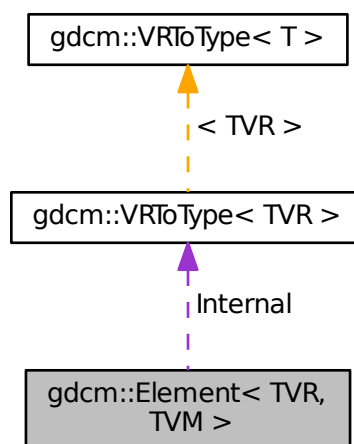

## Public Types

- typedef [VRToType](#)< TVR >::Type Type

## Public Member Functions

- [DataElement](#) [GetAsDataElement](#) () const
- unsigned long [GetLength](#) () const
- const [VRToType](#)< TVR >::Type & [GetValue](#) (unsigned int idx=0) const
- [VRToType](#)< TVR >::Type & [GetValue](#) (unsigned int idx=0)
- const [VRToType](#)< TVR >::Type \* [GetValues](#) () const
- [VRToType](#)< TVR >::Type [operator\[\]](#) (unsigned int idx) const
- void [Print](#) (std::ostream &\_os) const
- void [Read](#) (std::istream &\_is)
- void [Set](#) ([Value](#) const &v)
- void [SetFromDataElement](#) ([DataElement](#) const &de)
- void [SetValue](#) (typename [VRToType](#)< TVR >::Type v, unsigned int idx=0)
- void [Write](#) (std::ostream &\_os) const

## Static Public Member Functions

- static [VM](#) [GetVM](#) ()
- static [VR](#) [GetVR](#) ()

## Public Attributes

- [VRToType](#)< TVR >::Type [Internal](#) [[VMToLength](#)< TVM >::Length]

## Protected Member Functions

- void [SetNoSwap](#) ([Value](#) const &v)

### 27.93.1 Detailed Description

template<int TVR, int TVM>class gdcmm::Element< TVR, TVM >

[Element](#) class.

Note

TODO

Examples:

[csa2img.cxx](#), [DumpADAC.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GetSubSequenceData.cxx](#), and [iU22tomultisc.cxx](#).

## 27.93.2 Member Typedef Documentation

27.93.2.1 `template<int TVR, int TVM> typedef VRToType<TVR>::Type gdcmm::Element< TVR, TVM >::Type`

## 27.93.3 Member Function Documentation

27.93.3.1 `template<int TVR, int TVM> DataElement gdcmm::Element< TVR, TVM >::GetAsDataElement ( ) const`  
[inline]

27.93.3.2 `template<int TVR, int TVM> unsigned long gdcmm::Element< TVR, TVM >::GetLength ( ) const` [inline]

27.93.3.3 `template<int TVR, int TVM> const VRToType<TVR>::Type& gdcmm::Element< TVR, TVM >::GetValue ( unsigned int idx = 0 ) const` [inline]

27.93.3.4 `template<int TVR, int TVM> VRToType<TVR>::Type& gdcmm::Element< TVR, TVM >::GetValue ( unsigned int idx = 0 )` [inline]

27.93.3.5 `template<int TVR, int TVM> const VRToType<TVR>::Type* gdcmm::Element< TVR, TVM >::GetValues ( ) const`  
[inline]

27.93.3.6 `template<int TVR, int TVM> static VM gdcmm::Element< TVR, TVM >::GetVM ( )` [inline],[static]

27.93.3.7 `template<int TVR, int TVM> static VR gdcmm::Element< TVR, TVM >::GetVR ( )` [inline],[static]

27.93.3.8 `template<int TVR, int TVM> VRToType<TVR>::Type gdcmm::Element< TVR, TVM >::operator[] ( unsigned int idx ) const` [inline]

27.93.3.9 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::Print ( std::ostream &_os ) const` [inline]

27.93.3.10 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::Read ( std::istream &_is )` [inline]

27.93.3.11 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::Set ( Value const & v )` [inline]

27.93.3.12 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::SetFromDataElement ( DataElement< TVR, TVM > const & de )` [inline]

27.93.3.13 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::SetNoSwap ( Value const & v )` [inline],  
[protected]

27.93.3.14 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::SetValue ( typename VRToType< TVR >::Type v, unsigned int idx = 0 )` [inline]

27.93.3.15 `template<int TVR, int TVM> void gdcmm::Element< TVR, TVM >::Write ( std::ostream &_os ) const` [inline]

## 27.93.4 Member Data Documentation

27.93.4.1 `template<int TVR, int TVM> VRToType<TVR>::Type gdcmm::Element< TVR, TVM >::Internal[VMToLength< TVM >::Length]`


The documentation for this class was generated from the following file:

- [gdcmmElement.h](#)


## 27.94 gdcmm::Element< TVR, VM::VM1\_2 > Class Template Reference

```
#include <gdcmmElement.h>
```

Inheritance diagram for gdcmm::Element< TVR, VM::VM1\_2 >:

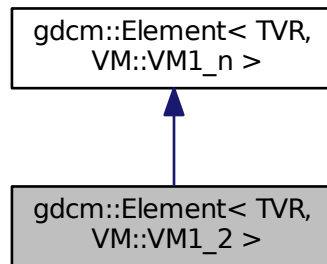

Collaboration diagram for gdcmm::Element< TVR, VM::VM1\_2 >:

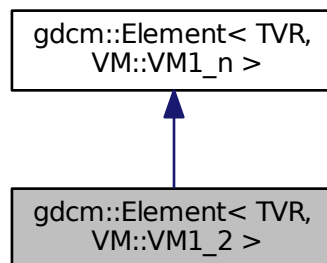

### Public Types

- typedef `Element< TVR, VM::VM1_n >` `Parent`

### Public Member Functions

- void `SetLength` (int len)

### Additional Inherited Members

### 27.94.1 Member Typedef Documentation

27.94.1.1 `template<int TVR> typedef Element<TVR, VM::VM1_n> gdcm::Element< TVR, VM::VM1_2 >::Parent`

### 27.94.2 Member Function Documentation

27.94.2.1 `template<int TVR> void gdcm::Element< TVR, VM::VM1_2 >::SetLength( int len ) [inline]`


The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.95 gdcm::Element< TVR, VM::VM1\_n > Class Template Reference

```
#include <gdcmElement.h>
```

Inheritance diagram for `gdcm::Element< TVR, VM::VM1_n >`:

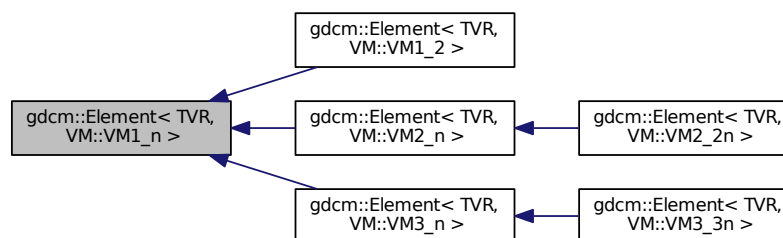

### Public Types

- typedef `VRToType< TVR >::Type` `Type`

### Public Member Functions

- `Element()`
- `Element(const Element &_val)`
- `~Element()`
- `DataElement GetAsDataElement() const`
- `unsigned long GetLength() const`
- `const VRToType< TVR >::Type & GetValue(unsigned int idx=0) const`
- `VRToType< TVR >::Type & GetValue(unsigned int idx=0)`
- `Element & operator=(const Element &_val)`
- `VRToType< TVR >::Type operator[](unsigned int idx) const`
- `void Print(std::ostream &_os) const`
- `void Read(std::istream &_is)`
- `void Set(Value const &v)`

- void [SetArray](#) (const [Type](#) \*array, unsigned long len, bool save=false)
- void [SetFromDataElement](#) ([DataElement](#) const &de)
- void [SetLength](#) (unsigned long len)
- void [SetValue](#) (typename [VRToType](#)< TVR >::Type v, unsigned int idx=0)
- void [Write](#) (std::ostream &\_os) const
- void [WriteASCII](#) (std::ostream &os) const

## Static Public Member Functions

- static [VM GetVM](#) ()
- static [VR GetVR](#) ()

## Protected Member Functions

- void [SetNoSwap](#) ([Value](#) const &v)

## 27.95.1 Member Typedef Documentation

27.95.1.1 `template<int TVR> typedef VRToType<TVR>::Type gdcm::Element< TVR, VM::VM1_n >::Type`

## 27.95.2 Constructor & Destructor Documentation

27.95.2.1 `template<int TVR> gdcm::Element< TVR, VM::VM1_n >::Element ( ) [inline], [explicit]`

27.95.2.2 `template<int TVR> gdcm::Element< TVR, VM::VM1_n >::~~Element ( ) [inline]`

27.95.2.3 `template<int TVR> gdcm::Element< TVR, VM::VM1_n >::Element ( const Element< TVR, VM::VM1_n > &_val ) [inline]`

## 27.95.3 Member Function Documentation

27.95.3.1 `template<int TVR> DataElement gdcm::Element< TVR, VM::VM1_n >::GetAsDataElement ( ) const [inline]`

References `gdcm::DataElement::GetVR()`, `gdcm::DataElement::SetByteValue()`, `gdcm::DataElement::SetVR()`, `gdcm::VR::SQ`, `gdcm::VR::UI`, and `gdcm::VR::VRASCII`.

27.95.3.2 `template<int TVR> unsigned long gdcm::Element< TVR, VM::VM1_n >::GetLength ( ) const [inline]`

27.95.3.3 `template<int TVR> const VRToType<TVR>::Type& gdcm::Element< TVR, VM::VM1_n >::GetValue ( unsigned int idx = 0 ) const [inline]`

27.95.3.4 `template<int TVR> VRToType<TVR>::Type& gdcm::Element< TVR, VM::VM1_n >::GetValue ( unsigned int idx = 0 ) [inline]`

27.95.3.5 `template<int TVR> static VM gdcm::Element< TVR, VM::VM1_n >::GetVM ( ) [inline], [static]`

References `gdcm::VM::VM1_n`.



27.95.3.6 `template<int TVR> static VR gdcm::Element< TVR, VM::VM1_n >::GetVR ( ) [inline],[static]`

27.95.3.7 `template<int TVR> Element& gdcm::Element< TVR, VM::VM1_n >::operator= ( const Element< TVR, VM::VM1_n > &_val ) [inline]`

27.95.3.8 `template<int TVR> VRToType<TVR>::Type gdcm::Element< TVR, VM::VM1_n >::operator[] ( unsigned int idx ) const [inline]`

27.95.3.9 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::Print ( std::ostream &_os ) const [inline]`

27.95.3.10 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::Read ( std::istream &_is ) [inline]`

27.95.3.11 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::Set ( Value const & v ) [inline]`

References `gdcm::ByteValue::GetLength()`, `gdcm::ByteValue::GetPointer()`, and `gdcm::VR::VRBINARY`.

27.95.3.12 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::SetArray ( const Type * array, unsigned long len, bool save = false ) [inline]`

27.95.3.13 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::SetFromDataElement ( DataElement< TVR, VM::VM1_n > const & de ) [inline]`

References `gdcm::DataElement::GetByteValue()`, `gdcm::DataElement::GetValue()`, `gdcm::DataElement::GetVR()`, `gdcm::VR::INVALID`, and `gdcm::VR::UN`.

27.95.3.14 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::SetLength ( unsigned long len ) [inline]`

27.95.3.15 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::SetNoSwap ( Value const & v ) [inline],[protected]`

References `gdcm::ByteValue::GetLength()`, `gdcm::ByteValue::GetPointer()`, and `gdcm::VR::VRBINARY`.

27.95.3.16 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::SetValue ( typename VRToType< TVR >::Type v, unsigned int idx = 0 ) [inline]`

27.95.3.17 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::Write ( std::ostream &_os ) const [inline]`

27.95.3.18 `template<int TVR> void gdcm::Element< TVR, VM::VM1_n >::WriteASCII ( std::ostream & os ) const [inline]`


The documentation for this class was generated from the following file:

- [gdcmElement.h](#)


## 27.96 gdcm::Element< TVR, VM::VM2\_2n > Class Template Reference

```
#include <gdcmElement.h>
```

Inheritance diagram for `gdcm::Element< TVR, VM::VM2_2n >`:

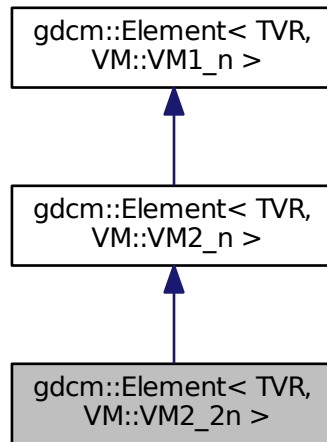

Collaboration diagram for `gdcm::Element< TVR, VM::VM2_2n >`:

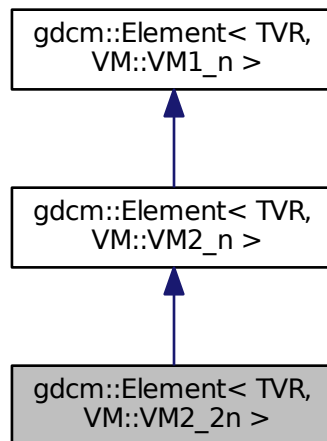

## Public Types

- typedef `Element< TVR, VM::VM2_n >` `Parent`

## Public Member Functions

- void [SetLength](#) (int len)

## Additional Inherited Members

### 27.96.1 Member Typedef Documentation

27.96.1.1 `template<int TVR> typedef Element<TVR, VM::VM2_n> gdcm::Element< TVR, VM::VM2_2n >::Parent`

### 27.96.2 Member Function Documentation

27.96.2.1 `template<int TVR> void gdcm::Element< TVR, VM::VM2_2n >::SetLength ( int len ) [inline]`


The documentation for this class was generated from the following file:

- [gdcmElement.h](#)


## 27.97 `gdcm::Element< TVR, VM::VM2_n >` Class Template Reference

```
#include <gdcmElement.h>
```

Inheritance diagram for `gdcm::Element< TVR, VM::VM2_n >`:

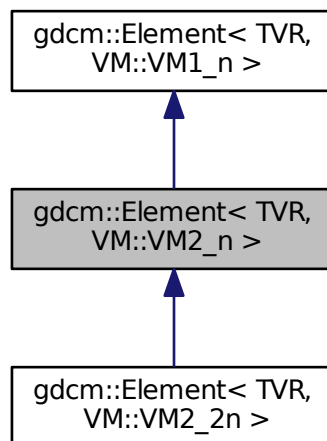

Collaboration diagram for `gdcM::Element< TVR, VM::VM2_n >`:

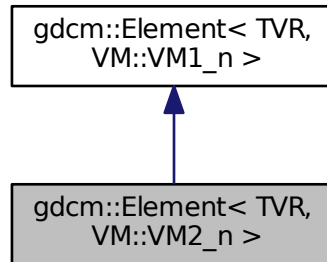

## Public Types

- typedef `Element< TVR, VM::VM1_n >` `Parent`

## Public Member Functions

- void `SetLength` (int *len*)

## Additional Inherited Members

### 27.97.1 Member Typedef Documentation

27.97.1.1 `template<int TVR> typedef Element<TVR, VM::VM1_n> gdcM::Element< TVR, VM::VM2_n >::Parent`

### 27.97.2 Member Function Documentation

27.97.2.1 `template<int TVR> void gdcM::Element< TVR, VM::VM2_n >::SetLength ( int len ) [inline]`


The documentation for this class was generated from the following file:

- `gdcMElement.h`


## 27.98 `gdcM::Element< TVR, VM::VM3_3n >` Class Template Reference

```
#include <gdcMElement.h>
```

Inheritance diagram for gdcm::Element< TVR, VM::VM3\_3n >:

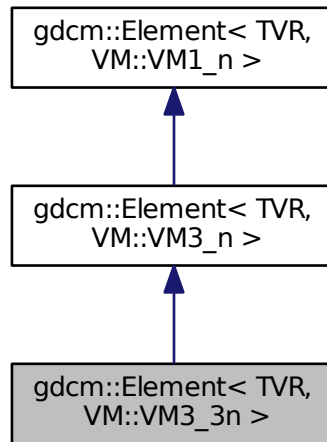

Collaboration diagram for gdcm::Element< TVR, VM::VM3\_3n >:

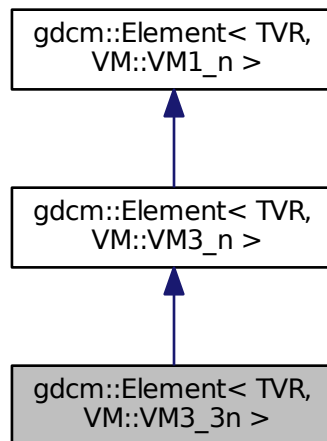

## Public Types

- typedef [Element](#)< TVR, [VM::VM3\\_n](#) > [Parent](#)

## Public Member Functions

- void [SetLength](#) (int len)

## Additional Inherited Members

### 27.98.1 Member Typedef Documentation

27.98.1.1 `template<int TVR> typedef Element<TVR, VM::VM3_n> gdcmm::Element< TVR, VM::VM3_3n >::Parent`

### 27.98.2 Member Function Documentation

27.98.2.1 `template<int TVR> void gdcmm::Element< TVR, VM::VM3_3n >::SetLength ( int len ) [inline]`


The documentation for this class was generated from the following file:

- [gdcmmElement.h](#)


## 27.99 gdcmm::Element< TVR, VM::VM3\_n > Class Template Reference

```
#include <gdcmmElement.h>
```

Inheritance diagram for gdcmm::Element< TVR, VM::VM3\_n >:

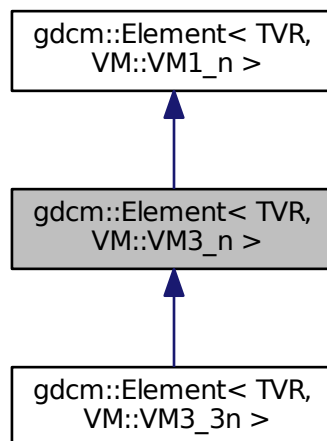

Collaboration diagram for `gdcm::Element< TVR, VM::VM3_n >`:

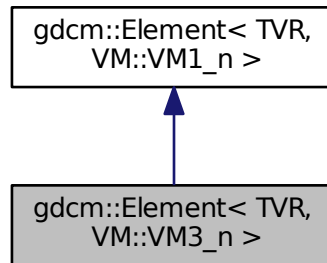

## Public Types

- typedef `Element< TVR, VM::VM1_n >` `Parent`

## Public Member Functions

- void `SetLength` (int len)

## Additional Inherited Members

### 27.99.1 Member Typedef Documentation

27.99.1.1 `template<int TVR> typedef Element<TVR, VM::VM1_n> gdcm::Element< TVR, VM::VM3_n >::Parent`

### 27.99.2 Member Function Documentation

27.99.2.1 `template<int TVR> void gdcm::Element< TVR, VM::VM3_n >::SetLength ( int len ) [inline]`

The documentation for this class was generated from the following file:

- `gdcmElement.h`

## 27.100 `gdcm::Element< VR::AS, VM::VM5 >` Class Template Reference

```
#include <gdcmElement.h>
```

## Public Member Functions

- unsigned long `GetLength` () const
- void `Print` (std::ostream &\_os) const

## Public Attributes

- char [Internal](#) [[VMToLength](#)< [VM::VM5](#) >::Length \*sizeof([VRToType](#)< [VR::AS](#) >::Type)]

## 27.100.1 Member Function Documentation

27.100.1.1 unsigned long [gdcm::Element](#)< [VR::AS](#), [VM::VM5](#) >::GetLength ( ) const [\[inline\]](#)

27.100.1.2 void [gdcm::Element](#)< [VR::AS](#), [VM::VM5](#) >::Print ( std::ostream &\_os ) const [\[inline\]](#)

## 27.100.2 Member Data Documentation

27.100.2.1 char [gdcm::Element](#)< [VR::AS](#), [VM::VM5](#) >::Internal[[VMToLength](#)< [VM::VM5](#) >::Length \*sizeof([VRToType](#)< [VR::AS](#) >::Type)]


The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.101 [gdcm::Element](#)< [VR::OB](#), [VM::VM1](#) > Class Template Reference


```
#include <gdcmElement.h>
```

Inheritance diagram for [gdcm::Element](#)< [VR::OB](#), [VM::VM1](#) >:

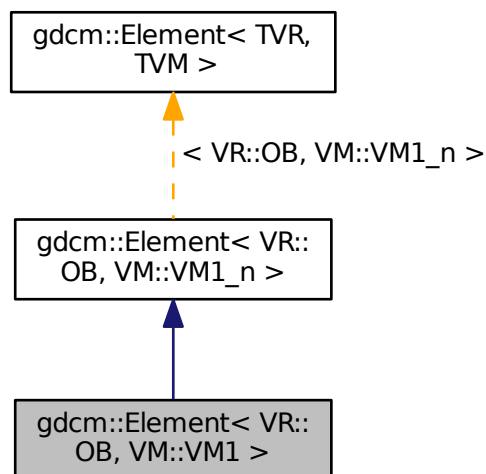



Collaboration diagram for gdcM::Element< VR::OB, VM::VM1 >:

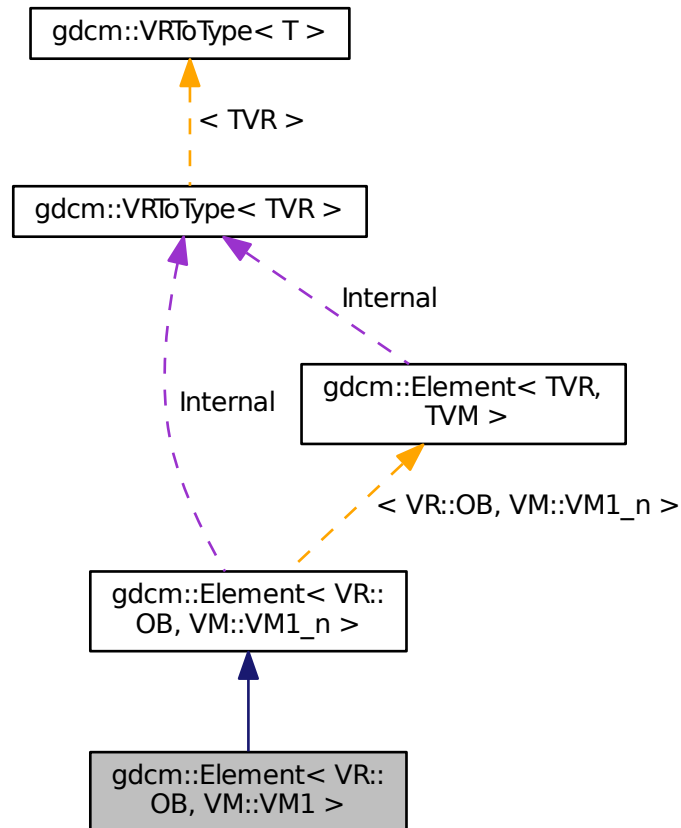

### Additional Inherited Members


The documentation for this class was generated from the following file:

- [gdcMElement.h](#)


## 27.102 gdcM::Element< VR::OW, VM::VM1 > Class Template Reference

```
#include <gdcMElement.h>
```

Inheritance diagram for gdcM::Element< VR::OW, VM::VM1 >:

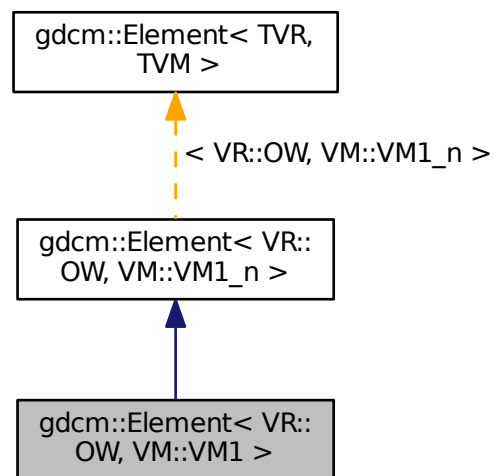

Collaboration diagram for gdcm::Element< VR::OW, VM::VM1 >:

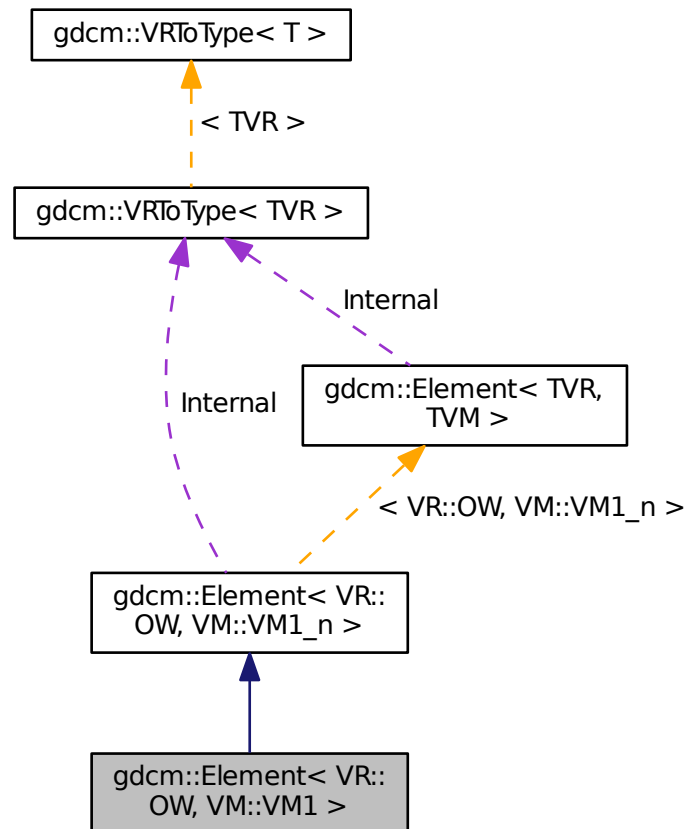

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.103 gdcm::ElementDisableCombinations< TVR, TVM > Class Template Reference

A class which is used to produce compile errors for an invalid combination of template parameters.

```
#include <gdcmElement.h>
```

### 27.103.1 Detailed Description

```
template<int TVR, int TVM>class gdcm::ElementDisableCombinations< TVR, TVM >
```

A class which is used to produce compile errors for an invalid combination of template parameters.

Invalid combinations have specialized declarations with no definition.

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.104 `gdcm::ElementDisableCombinations< VR::OB, VM::VM1_n >` Class Template Reference

```
#include <gdcmElement.h>
```

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.105 `gdcm::ElementDisableCombinations< VR::OW, VM::VM1_n >` Class Template Reference

```
#include <gdcmElement.h>
```

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.106 `gdcm::EncapsulatedDocument` Class Reference

[EncapsulatedDocument.](#)

```
#include <gdcmEncapsulatedDocument.h>
```

### Public Member Functions

- [EncapsulatedDocument\(\)](#)

### 27.106.1 Detailed Description

[EncapsulatedDocument.](#)

### 27.106.2 Constructor & Destructor Documentation

27.106.2.1 `gdcm::EncapsulatedDocument::EncapsulatedDocument ( )` `[inline]`

The documentation for this class was generated from the following file:

- [gdcmEncapsulatedDocument.h](#)

## 27.107 gdcm::EncodingImplementation< T > Class Template Reference

[EncodingImplementation](#).

```
#include <gdcmElement.h>
```

### 27.107.1 Detailed Description

```
template<int T>class gdcm::EncodingImplementation< T >
```

[EncodingImplementation](#).

Note

TODO

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.108 gdcm::EncodingImplementation< VR::VRASCII > Class Template Reference

```
#include <gdcmElement.h>
```

### Public Member Functions

- template<>  
void [Write](#) (const float \*data, unsigned long length, std::ostream &\_os)
- template<>  
void [Write](#) (const double \*data, unsigned long length, std::ostream &\_os)

### Static Public Member Functions

- template<typename T >  
static void [Read](#) (T \*data, unsigned long length, std::istream &\_is)
- template<typename T >  
static void [ReadComputeLength](#) (T \*data, unsigned int &length, std::istream &\_is)
- template<typename T >  
static void [ReadNoSwap](#) (T \*data, unsigned long length, std::istream &\_is)
- template<typename T >  
static void [Write](#) (const T \*data, unsigned long length, std::ostream &\_os)

### 27.108.1 Member Function Documentation

27.108.1.1 template<typename T > static void gdcm::EncodingImplementation< VR::VRASCII >::Read ( T \* data, unsigned long length, std::istream &\_is ) [inline], [static]

27.108.1.2 `template<typename T> static void gdcm::EncodingImplementation< VR::VRASCII >::ReadComputeLength ( T * data, unsigned int & length, std::istream & _is ) [inline], [static]`

References `gdcm::backslash()`.

27.108.1.3 `template<typename T> static void gdcm::EncodingImplementation< VR::VRASCII >::ReadNoSwap ( T * data, unsigned long length, std::istream & _is ) [inline], [static]`

27.108.1.4 `template<typename T> static void gdcm::EncodingImplementation< VR::VRASCII >::Write ( const T * data, unsigned long length, std::ostream & _os ) [inline], [static]`

27.108.1.5 `template<> void gdcm::EncodingImplementation< VR::VRASCII >::Write ( const float * data, unsigned long length, std::ostream & _os ) [inline]`

References `gdcm::to_string()`.

27.108.1.6 `template<> void gdcm::EncodingImplementation< VR::VRASCII >::Write ( const double * data, unsigned long length, std::ostream & _os ) [inline]`

References `gdcm::to_string()`.

The documentation for this class was generated from the following file:

- [gdcmElement.h](#)

## 27.109 gdcm::EncodingImplementation< VR::VRBINARY > Class Template Reference

```
#include <gdcmElement.h>
```

### Static Public Member Functions

- `template<typename T> static void Read (T *data, unsigned long length, std::istream &_is)`
- `template<typename T> static void ReadComputeLength (T *data, unsigned int &length, std::istream &_is)`
- `template<typename T> static void ReadNoSwap (T *data, unsigned long length, std::istream &_is)`
- `template<typename T> static void Write (const T *data, unsigned long length, std::ostream &_os)`

### 27.109.1 Member Function Documentation

27.109.1.1 `template<typename T> static void gdcm::EncodingImplementation< VR::VRBINARY >::Read ( T * data, unsigned long length, std::istream & _is ) [inline], [static]`

References `gdcm::SwapperNoOp::SwapArray()`.

27.109.1.2 `template<typename T> static void gdcm::EncodingImplementation< VR::VRBINARY >::ReadComputeLength ( T * data, unsigned int & length, std::istream & _is ) [inline], [static]`

27.109.1.3 `template<typename T> static void gdcm::EncodingImplementation< VR::VRBINARY >::ReadNoSwap ( T * data, unsigned long length, std::istream & _is ) [inline], [static]`

27.109.1.4 `template<typename T> static void gdcm::EncodingImplementation< VR::VRBINARY >::Write ( const T * data, unsigned long length, std::ostream & _os ) [inline], [static]`

References `gdcm::SwapperNoOp::Swap()`.


The documentation for this class was generated from the following file:

- [gdcmElement.h](#)


## 27.110 gdcm::EndEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::EndEvent`:

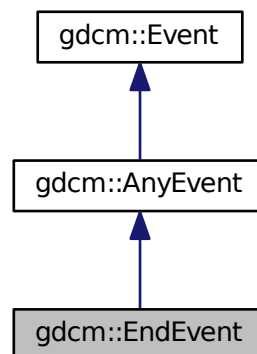

Collaboration diagram for `gdcm::EndEvent`:

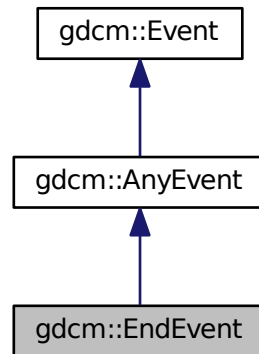

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.111 gdcm::EnumeratedValues Class Reference

**Element.** A Data [Element](#) with Enumerated Values that does not have a [Value](#) equivalent to one of the Values specified in this standard has an invalid value within the scope of a specific Information Object/SOP Class definition. Note:

```
#include <gdcmEnumeratedValues.h>
```

### Public Member Functions

- [EnumeratedValues](#) ()

#### 27.111.1 Detailed Description

**Element.** A Data [Element](#) with Enumerated Values that does not have a [Value](#) equivalent to one of the Values specified in this standard has an invalid value within the scope of a specific Information Object/SOP Class definition. Note:

1. [Patient](#) Sex (0010, 0040) is an example of a Data [Element](#) having Enumerated Values. It is defined to have a [Value](#) that is either "M", "F", or "O" (see PS 3.3). No other [Value](#) shall be given to this Data [Element](#).
2. Future modifications of this standard may add to the set of allowed values for Data Elements with Enumerated Values. Such additions by themselves may or may not require a change in SOP Class [UIDs](#), depending on the semantics of the Data [Element](#).



### 27.111.2 Constructor & Destructor Documentation

#### 27.111.2.1 gdcm::EnumeratedValues::EnumeratedValues ( ) [inline]

The documentation for this class was generated from the following file:


- [gdcmEnumeratedValues.h](#)

## 27.112 gdcm::Event Class Reference

superclass for callback/observer methods

```
#include <gdcmEvent.h>
```

Inheritance diagram for gdcm::Event:

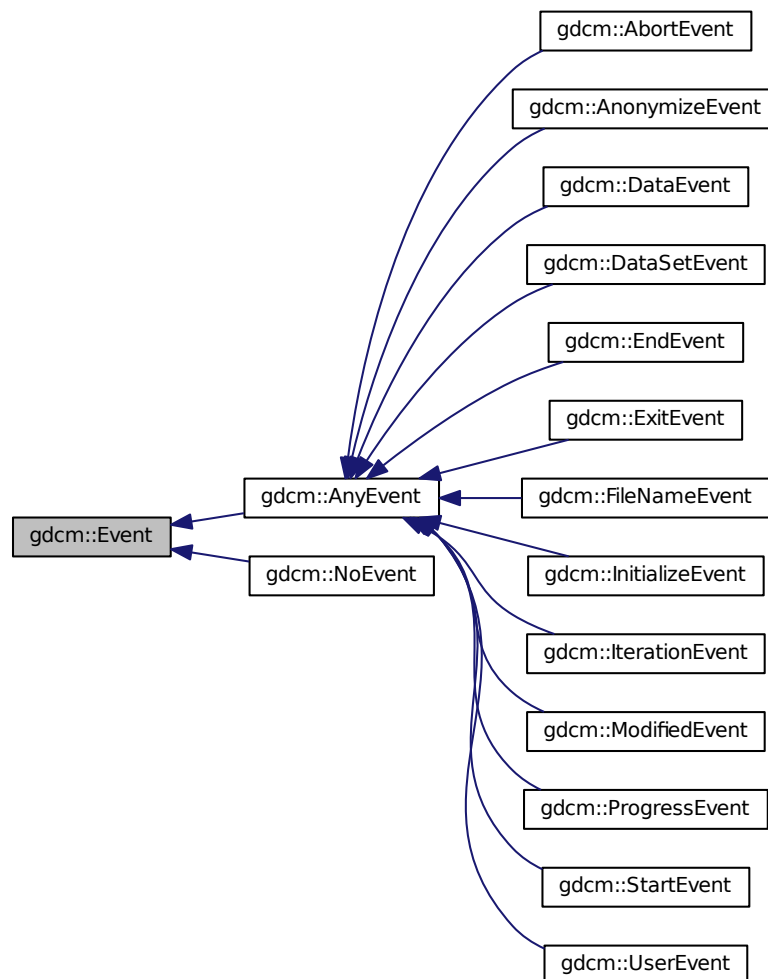

## Public Member Functions

- [Event](#) ()
- [Event](#) (const [Event](#) &)
- virtual [~Event](#) ()
- virtual bool [CheckEvent](#) (const [Event](#) \*) const =0
- virtual const char \* [GetEventName](#) (void) const =0
- virtual [Event](#) \* [MakeObject](#) () const =0
- virtual void [Print](#) (std::ostream &os) const

### 27.112.1 Detailed Description

superclass for callback/observer methods

See also

[Command Subject](#)

Examples:

[SimpleScanner.cxx](#).

### 27.112.2 Constructor & Destructor Documentation

27.112.2.1 `gdcm::Event::Event ( )`

27.112.2.2 `gdcm::Event::Event ( const Event & )`

27.112.2.3 `virtual gdcm::Event::~~Event ( )` `[virtual]`

### 27.112.3 Member Function Documentation

27.112.3.1 `virtual bool gdcm::Event::CheckEvent ( const Event * ) const` `[pure virtual]`

Check if given event matches or derives from this event.

27.112.3.2 `virtual const char* gdcm::Event::GetEventName ( void ) const` `[pure virtual]`

Return the StringName associated with the event.

Implemented in [gdcm::FileNameEvent](#), [gdcm::ProgressEvent](#), [gdcm::DataSetEvent](#), [gdcm::AnonymizeEvent](#), and [gdcm::DataEvent](#).

27.112.3.3 `virtual Event* gdcm::Event::MakeObject ( ) const` `[pure virtual]`

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implemented in [gdcm::FileNameEvent](#), [gdcm::ProgressEvent](#), [gdcm::DataSetEvent](#), [gdcm::AnonymizeEvent](#), and [gdcm::DataEvent](#).

27.112.3.4 virtual void gdcM::Event::Print ( std::ostream & os ) const [virtual]

Print [Event](#) information. This method can be overridden by specific [Event](#) subtypes. The default is to print out the type of the event.

Referenced by gdcM::operator<<().

The documentation for this class was generated from the following file:


- [gdcMEvent.h](#)

## 27.113 gdcM::Exception Class Reference


[Exception](#).

```
#include <gdcMException.h>
```

Inheritance diagram for gdcM::Exception:

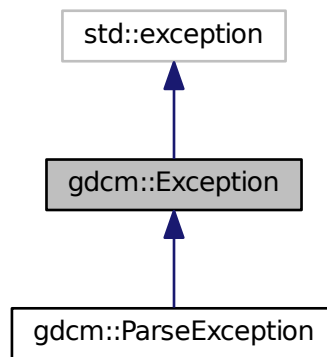

Collaboration diagram for gdcM::Exception:

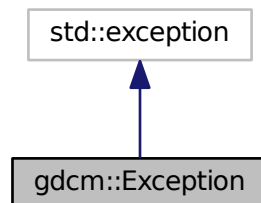

## Public Member Functions

- [Exception](#) (const char \*desc="None", const char \*file=\_\_FILE\_\_, unsigned int lineNumber=\_\_LINE\_\_, const char \*func="")
- virtual [~Exception](#) () throw ()
- const char \* [GetDescription](#) () const  
*Return the Description.*
- const char \* [what](#) () const throw ()  
*what implementation*

### 27.113.1 Detailed Description

[Exception](#).

Standard exception handling object.

#### Note

Its copy-constructor and assignment operator are generated by the compiler.

### 27.113.2 Constructor & Destructor Documentation

**27.113.2.1** `gdcmm::Exception::Exception ( const char * desc = "None", const char * file = __FILE__, unsigned int lineNumber = __LINE__, const char * func = " " ) [inline],[explicit]`

Explicit constructor, initializing the description and the text returned by [what\(\)](#).

#### Note

The last parameter is ignored for the time being. It may be used to specify the function where the exception was thrown.

**27.113.2.2** `virtual gdcmm::Exception::~~Exception ( ) throw ) [inline],[virtual]`

### 27.113.3 Member Function Documentation

**27.113.3.1** `const char* gdcmm::Exception::GetDescription ( ) const [inline]`

Return the Description.

Referenced by `gdcmm::SequenceOfItems::Read()`.

**27.113.3.2** `const char* gdcmm::Exception::what ( ) const throw ) [inline]`

*what implementation*

Referenced by `gdcmm::SequenceOfFragments::ReadValue()`.


The documentation for this class was generated from the following file:

- [gdcmmException.h](#)


## 27.114 gdcm::ExitEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for gdcm::ExitEvent:

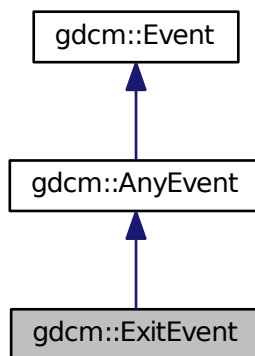

Collaboration diagram for gdcm::ExitEvent:

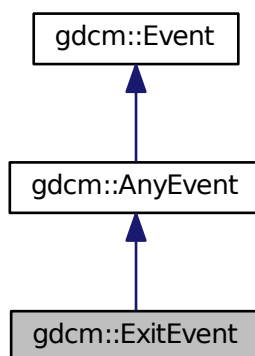

### Additional Inherited Members

The documentation for this class was generated from the following file:


- [gdcmEvent.h](#)

## 27.115 gdcm::ExplicitDataElement Class Reference


Class to read/write a [DataElement](#) as Explicit Data [Element](#).

```
#include <gdcmExplicitDataElement.h>
```

Inheritance diagram for gdcm::ExplicitDataElement:

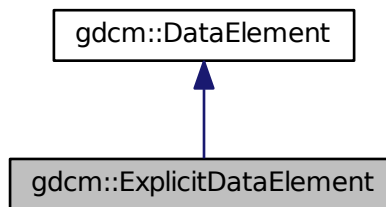

Collaboration diagram for gdcm::ExplicitDataElement:

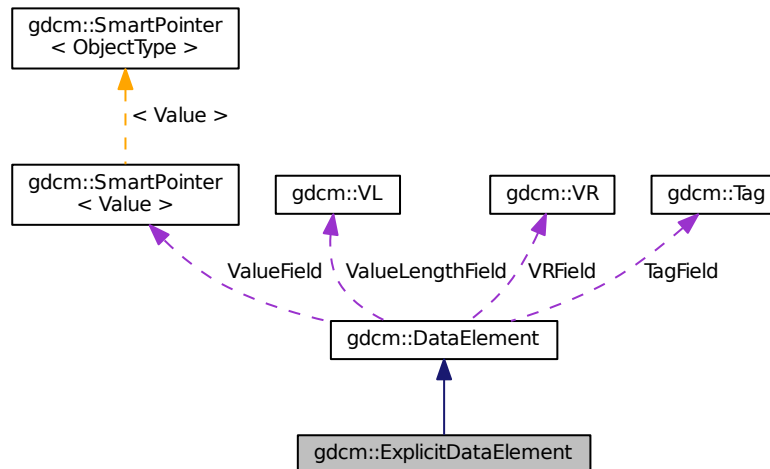

### Public Member Functions

- [VL GetLength](#) () const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadPreValue](#) (std::istream &is)

- `template<typename TSwap >`  
`std::istream & ReadValue (std::istream &is, bool readvalues=true)`
- `template<typename TSwap >`  
`std::istream & ReadWithLength (std::istream &is, VL &length)`
- `template<typename TSwap >`  
`const std::ostream & Write (std::ostream &os) const`

## Additional Inherited Members

### 27.115.1 Detailed Description

Class to read/write a [DataElement](#) as Explicit Data [Element](#).

#### Note

bla

### 27.115.2 Member Function Documentation

27.115.2.1 `VL gdcm::ExplicitDataElement::GetLength ( ) const`

27.115.2.2 `template<typename TSwap > std::istream& gdcm::ExplicitDataElement::Read ( std::istream & is )`

27.115.2.3 `template<typename TSwap > std::istream& gdcm::ExplicitDataElement::ReadPreValue ( std::istream & is )`

27.115.2.4 `template<typename TSwap > std::istream& gdcm::ExplicitDataElement::ReadValue ( std::istream & is, bool readvalues = true )`

27.115.2.5 `template<typename TSwap > std::istream& gdcm::ExplicitDataElement::ReadWithLength ( std::istream & is, VL & length )`

27.115.2.6 `template<typename TSwap > const std::ostream& gdcm::ExplicitDataElement::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:


- [gdcmExplicitDataElement.h](#)

## 27.116 gdcm::ExplicitImplicitDataElement Class Reference


Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#).

```
#include <gdcmExplicitImplicitDataElement.h>
```

Inheritance diagram for `gdc::ExplicitImplicitDataElement`:

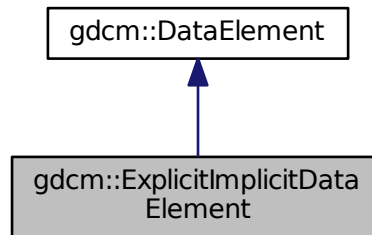

Collaboration diagram for `gdc::ExplicitImplicitDataElement`:

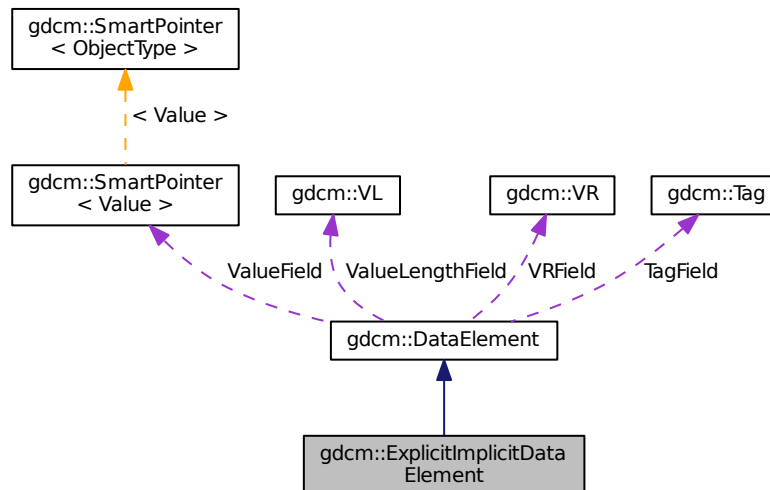

## Public Member Functions

- `VL GetLength () const`
- `template<typename TSwap > std::istream & Read (std::istream &is)`
- `template<typename TSwap > std::istream & ReadPreValue (std::istream &is)`
- `template<typename TSwap > std::istream & ReadValue (std::istream &is, bool readvalues=true)`
- `template<typename TSwap > std::istream & ReadWithLength (std::istream &is, VL &length)`



## Additional Inherited Members

### 27.116.1 Detailed Description

Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#).

#### Note

This only happen for some Philips images Should I derive from [ExplicitDataElement](#) instead ? This is the class that is the closest the GDCM1.x parser. At each element we try first to read it as explicit, if this fails, then we try again as an implicit element.

### 27.116.2 Member Function Documentation

27.116.2.1 `VL gdcm::ExplicitImplicitDataElement::GetLength ( ) const`

27.116.2.2 `template<typename TSwap > std::istream& gdcm::ExplicitImplicitDataElement::Read ( std::istream & is )`

27.116.2.3 `template<typename TSwap > std::istream& gdcm::ExplicitImplicitDataElement::ReadPreValue ( std::istream & is )`

27.116.2.4 `template<typename TSwap > std::istream& gdcm::ExplicitImplicitDataElement::ReadValue ( std::istream & is, bool readvalues = true )`

27.116.2.5 `template<typename TSwap > std::istream& gdcm::ExplicitImplicitDataElement::ReadWithLength ( std::istream & is, VL & length ) [inline]`

The documentation for this class was generated from the following file:

- [gdcmExplicitImplicitDataElement.h](#)

## 27.117 gdcm::Fiducials Class Reference

[Fiducials](#).

```
#include <gdcmFiducials.h>
```

### Public Member Functions

- [Fiducials](#) ()

### 27.117.1 Detailed Description

[Fiducials](#).

### 27.117.2 Constructor & Destructor Documentation

27.117.2.1 `gdcm::Fiducials::Fiducials ( ) [inline]`

The documentation for this class was generated from the following file:


- [gdcmFiducials.h](#)

## 27.118 gdcm::File Class Reference


a DICOM [File](#) See PS 3.10 [File](#): A [File](#) is an ordered string of zero or more bytes, where the first byte is at the beginning of the file and the last byte at the end of the [File](#). Files are identified by a unique [File](#) ID and may be written, read and/or deleted.

```
#include <gdcmFile.h>
```

Inheritance diagram for gdcm::File:

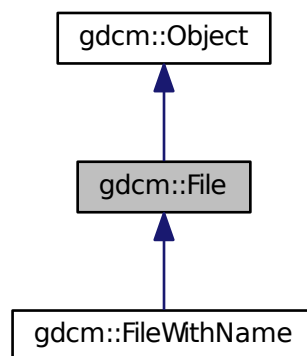

Collaboration diagram for gdcm::File:

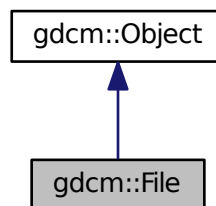

### Public Member Functions

- [File](#) ()
- const [DataSet](#) & [GetDataSet](#) () const

- Get Data Set.*
- [DataSet](#) & [GetDataSet](#) ()
- Get Data Set.*
- const [FileMetaInformation](#) & [GetHeader](#) () const
- Get File Meta Information.*
- [FileMetaInformation](#) & [GetHeader](#) ()
- Get File Meta Information.*
- std::istream & [Read](#) (std::istream &is)
- Read.*
- void [SetDataSet](#) (const [DataSet](#) &ds)
- Set Data Set.*
- void [SetHeader](#) (const [FileMetaInformation](#) &fmi)
- Set File Meta Information.*
- std::ostream const & [Write](#) (std::ostream &os) const
- Write.*

## Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [File](#) &val)

## Additional Inherited Members

### 27.118.1 Detailed Description

a DICOM [File](#) See PS 3.10 [File](#): A [File](#) is an ordered string of zero or more bytes, where the first byte is at the beginning of the file and the last byte at the end of the [File](#). Files are identified by a unique [File](#) ID and may be written, read and/or deleted.

See also

[Reader Writer](#)

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [DiffFile.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [EncapsulateFileInRawData.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenFakeImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [PatchFile.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadGEMSSDO.cxx](#), and [StreamImageReaderTest.cxx](#).

### 27.118.2 Constructor & Destructor Documentation

27.118.2.1 `gdcm::File::File ( ) [inline]`

### 27.118.3 Member Function Documentation

27.118.3.1 `const DataSet& gdcm::File::GetDataSet ( ) const [inline]`

Get Data Set.

**Examples:**

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [gdcmrtnplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [rle2img.cxx](#), and [StreamImageReaderTest.cxx](#).

**27.118.3.2 DataSet& gdcm::File::GetDataSet ( ) [inline]**

Get Data Set.

**27.118.3.3 const FileMetaInformation& gdcm::File::GetHeader ( ) const [inline]**

Get [File](#) Meta Information.

**Examples:**

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GetJPEGSamplePrecision.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [MpegVideoInfo.cs](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [rle2img.cxx](#), and [StreamImageReaderTest.cxx](#).

Referenced by `gdcm::operator<<()`.

**27.118.3.4 FileMetaInformation& gdcm::File::GetHeader ( ) [inline]**

Get [File](#) Meta Information.

**27.118.3.5 std::istream& gdcm::File::Read ( std::istream & is )**

Read.

**27.118.3.6 void gdcm::File::SetDataSet ( const DataSet & ds ) [inline]**

Set Data Set.

**27.118.3.7 void gdcm::File::SetHeader ( const FileMetaInformation & fmi ) [inline]**

Set [File](#) Meta Information.

**27.118.3.8 std::ostream const& gdcm::File::Write ( std::ostream & os ) const**

Write.

## 27.118.4 Friends And Related Function Documentation

27.118.4.1 `std::ostream& operator<< ( std::ostream & os, const File & val )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmFile.h](#)

## 27.119 gdcm::FileAnonymizer Class Reference


[FileAnonymizer](#).

```
#include <gdcmFileAnonymizer.h>
```

Inheritance diagram for gdcm::FileAnonymizer:

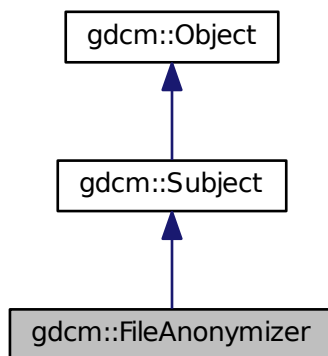

Collaboration diagram for `gdcm::FileAnonymizer`:

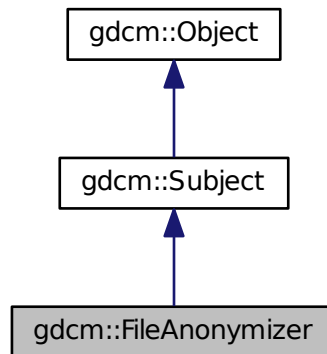

## Public Member Functions

- [FileAnonymizer](#) ()
- [~FileAnonymizer](#) ()
- void [Empty](#) ([Tag](#) const &t)
- void [Remove](#) ([Tag](#) const &t)  
*remove a tag (even a SQ can be removed)*
- void [Replace](#) ([Tag](#) const &t, const char \*value\_str)
- void [Replace](#) ([Tag](#) const &t, const char \*value\_data, [VL](#) const &vl)
- void [SetInputFileName](#) (const char \*filename\_native)  
*Set input filename.*
- void [SetOutputFileName](#) (const char \*filename\_native)  
*Set output filename.*
- bool [Write](#) ()  
*Write the output file.*

## Additional Inherited Members

### 27.119.1 Detailed Description

[FileAnonymizer](#).

This [Anonymizer](#) is a file-based [Anonymizer](#). It requires a valid DICOM file and will use the [Value](#) Length to skip over any information.

It will not load the DICOM dataset taken from [SetInputFileName\(\)](#) into memory and should consume much less memory than [Anonymizer](#).

## Warning

: Each time you call [Replace\(\)](#) with a value. This value will be copied, and stored in memory. The behavior is not ideal for extremely large data (larger than memory size). This class is really meant to take a large DICOM input file and then only change some small attribute.

## caveats:

- This class will NOT work with unordered attributes in a DICOM [File](#),
- This class does neither recompute nor update the Group Length element,
- This class currently does not update the [File](#) Meta Information header.
- Only strict inplace Replace operation is supported when input and output file are the same.

## 27.119.2 Constructor &amp; Destructor Documentation

27.119.2.1 `gdcm::FileAnonymizer::FileAnonymizer ( )`

27.119.2.2 `gdcm::FileAnonymizer::~~FileAnonymizer ( )`

## 27.119.3 Member Function Documentation

27.119.3.1 `void gdcm::FileAnonymizer::Empty ( Tag const & t )`

Make [Tag](#) t empty. Warning: does not handle SQ element

27.119.3.2 `void gdcm::FileAnonymizer::Remove ( Tag const & t )`

remove a tag (even a SQ can be removed)

27.119.3.3 `void gdcm::FileAnonymizer::Replace ( Tag const & t, const char * value_str )`

Replace tag with another value, if tag is not found it will be created: WARNING: this function can only execute if tag is a VRASCII. WARNING: Do not ever try to write a value in a SQ Data [Element](#) !

27.119.3.4 `void gdcm::FileAnonymizer::Replace ( Tag const & t, const char * value_data, VL const & vl )`

when the value contains \0, it is a good idea to specify the length. This function is required when dealing with VRBINARY tag

27.119.3.5 `void gdcm::FileAnonymizer::SetInputFileName ( const char * filename_native )`

Set input filename.

## Examples:

[FileAnonymize.cs](#).

27.119.3.6 void gdcmm::FileAnonymizer::SetOutputFileName ( const char \* *filename\_native* )

Set output filename.

27.119.3.7 bool gdcmm::FileAnonymizer::Write ( )

Write the output file.

The documentation for this class was generated from the following file:


- [gdcmmFileAnonymizer.h](#)

## 27.120 gdcmm::FileChangeTransferSyntax Class Reference

[FileChangeTransferSyntax](#).


```
#include <gdcmmFileChangeTransferSyntax.h>
```

Inheritance diagram for gdcmm::FileChangeTransferSyntax:

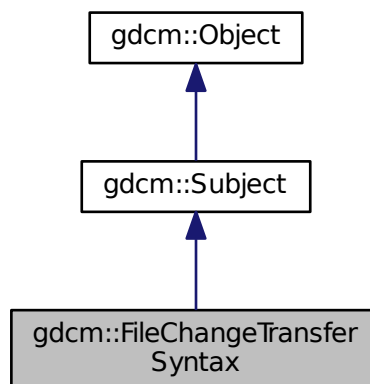



Collaboration diagram for gdcm::FileChangeTransferSyntax:

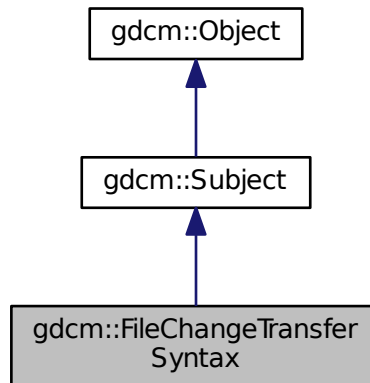

### Public Member Functions

- [FileChangeTransferSyntax](#) ()
- [~FileChangeTransferSyntax](#) ()
- [bool Change](#) ()  
*Change the transfer syntax.*
- [ImageCodec \\* GetCodec](#) ()
- [void SetInputFileName](#) (const char \*filename\_native)  
*Set input filename (raw DICOM)*
- [void SetOutputFileName](#) (const char \*filename\_native)  
*Set output filename (target compressed DICOM)*
- [void SetTransferSyntax](#) ([TransferSyntax](#) const &ts)  
*Specify the Target Transfer Syntax.*

### Static Public Member Functions

- static [SmartPointer< FileChangeTransferSyntax > New](#) ()  
*for wrapped language: instantiate a reference counted object*

### Additional Inherited Members

#### 27.120.1 Detailed Description

[FileChangeTransferSyntax](#).

This class is a file-based (limited) replacement of the in-memory [ImageChangeTransferSyntax](#).

This class provide a file-based compression-only mechanism. It will take in an uncompressed DICOM image file (Pixel Data element). Then produced as output a compressed DICOM file (Transfer Syntax will be updated).

Currently it supports the following transfer syntax:

- JPEGLosslessProcess14\_1

## 27.120.2 Constructor & Destructor Documentation

27.120.2.1 `gdcm::FileChangeTransferSyntax::FileChangeTransferSyntax ( )`

27.120.2.2 `gdcm::FileChangeTransferSyntax::~~FileChangeTransferSyntax ( )`

## 27.120.3 Member Function Documentation

27.120.3.1 `bool gdcm::FileChangeTransferSyntax::Change ( )`

Change the transfer syntax.

27.120.3.2 `ImageCodec* gdcm::FileChangeTransferSyntax::GetCodec ( )`

Retrieve the actual codec (valid after calling SetTransferSyntax) Only advanced users should call this function.

27.120.3.3 `static SmartPointer<FileChangeTransferSyntax> gdcm::FileChangeTransferSyntax::New ( ) [inline], [static]`

for wrapped language: instantiate a reference counted object

27.120.3.4 `void gdcm::FileChangeTransferSyntax::SetInputFileName ( const char * filename_native )`

Set input filename (raw DICOM)

27.120.3.5 `void gdcm::FileChangeTransferSyntax::SetOutputFileName ( const char * filename_native )`

Set output filename (target compressed DICOM)

27.120.3.6 `void gdcm::FileChangeTransferSyntax::SetTransferSyntax ( TransferSyntax const & ts )`

Specify the Target Transfer Syntax.

The documentation for this class was generated from the following file:

- [gdcmFileChangeTransferSyntax.h](#)

## 27.121 gdcm::FileDerivation Class Reference

[FileDerivation](#) class See PS 3.16 - 2008 For the list of Code [Value](#) that can be used for in Derivation Code Sequence.

```
#include <gdcmFileDerivation.h>
```

## Public Member Functions

- [FileDerivation](#) ()
- [~FileDerivation](#) ()
- bool [AddReference](#) (const char \*referencedsopclassuid, const char \*referencedsopinstanceuid)
- bool [Derive](#) ()
  - Change.*
- [File](#) & [GetFile](#) ()
- const [File](#) & [GetFile](#) () const
- void [SetDerivationCodeSequenceCodeValue](#) (unsigned int codevalue)
  - Specify the Derivation Code Sequence Code [Value](#). Eg 113040.*
- void [SetDerivationDescription](#) (const char \*dd)
  - Specify the Derivation Description. Eg "lossy conversion".*
- void [SetFile](#) (const [File](#) &f)
  - Set/Get [File](#).*
- void [SetPurposeOfReferenceCodeSequenceCodeValue](#) (unsigned int codevalue)
  - Specify the Purpose Of Reference Code [Value](#). Eg. 121320.*

## Protected Member Functions

- bool [AddDerivationDescription](#) ()
- bool [AddPurposeOfReferenceCodeSequence](#) ([DataSet](#) &ds)
- bool [AddSourceImageSequence](#) ()

### 27.121.1 Detailed Description

[FileDerivation](#) class See PS 3.16 - 2008 For the list of Code [Value](#) that can be used for in Derivation Code Sequence.

URL: [http://medical.nema.org/medical/dicom/2008/08\\_16pu.pdf](http://medical.nema.org/medical/dicom/2008/08_16pu.pdf)

DICOM Part 16 has two Context Groups CID 7202 and CID 7203 which contain a set of codes defining reason for a source image reference (ie. reason code for referenced image sequence) and a coded description of the deriation applied to the new image data from the original. Both these context groups are extensible.

[File](#) Derivation is compulsory when creating a lossy derived image.

Examples:

[GenFakelImage.cxx](#).

### 27.121.2 Constructor & Destructor Documentation

27.121.2.1 `gdcm::FileDerivation::FileDerivation ( )`

27.121.2.2 `gdcm::FileDerivation::~~FileDerivation ( )`

### 27.121.3 Member Function Documentation

27.121.3.1 `bool gdcm::FileDerivation::AddDerivationDescription ( )` [protected]

27.121.3.2 `bool gdcm::FileDerivation::AddPurposeOfReferenceCodeSequence ( DataSet & ds )` [protected]

27.121.3.3 `bool gdcM::FileDerivation::AddReference ( const char * referencedsopclassuid, const char * referencedsopinstanceuid )`

Create the proper reference. Need to pass the original SOP Class UID and the original SOP Instance UID, so that those value can be used as Reference.

Warning

`referencedsopclassuid` and `referencedsopinstanceuid` needs to be \0 padded. This is not compatible with how `ByteValue->GetPointer` works.

Examples:

[GenFakelImage.cxx](#).

27.121.3.4 `bool gdcM::FileDerivation::AddSourceImageSequence ( ) [protected]`

27.121.3.5 `bool gdcM::FileDerivation::Derive ( )`

Change.

Examples:

[GenFakelImage.cxx](#).

27.121.3.6 `File& gdcM::FileDerivation::GetFile ( ) [inline]`

Examples:

[GenFakelImage.cxx](#).

27.121.3.7 `const File& gdcM::FileDerivation::GetFile ( ) const [inline]`

27.121.3.8 `void gdcM::FileDerivation::SetDerivationCodeSequenceCodeValue ( unsigned int codevalue )`

Specify the Derivation Code Sequence Code [Value](#). Eg 113040.

Examples:

[GenFakelImage.cxx](#).

27.121.3.9 `void gdcM::FileDerivation::SetDerivationDescription ( const char * dd )`

Specify the Derivation Description. Eg "lossy conversion".

27.121.3.10 `void gdcM::FileDerivation::SetFile ( const File & f ) [inline]`

Set/Get [File](#).

Examples:

[GenFakelImage.cxx](#).

27.121.3.11 void gdcm::FileDerivation::SetPurposeOfReferenceCodeSequenceCodeValue ( unsigned int *codevalue* )

Specify the Purpose Of Reference Code [Value](#). Eg. 121320.

Examples:

[GenFakelImage.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmFileDerivation.h](#)

## 27.122 gdcm::FileExplicitFilter Class Reference

[FileExplicitFilter](#) class After changing a file from Implicit to Explicit representation (see [ImageChangeTransferSyntax](#)) one operation is to make sure the [VR](#) of each DICOM attribute are accurate and do match the one from PS 3.6. Indeed when a file is written in Implicit representation, the [VR](#) is not stored directly in the file.

```
#include <gdcmFileExplicitFilter.h>
```

### Public Member Functions

- [FileExplicitFilter](#) ()
- [~FileExplicitFilter](#) ()
- bool [Change](#) ()  
*Set FMI Transfer Syntax.*
- [File](#) & [GetFile](#) ()
- void [SetChangePrivateTags](#) (bool b)  
*Decide whether or not to [VR](#)ify private tags.*
- void [SetFile](#) (const [File](#) &f)  
*Set/Get [File](#).*
- void [SetRecomputeItemLength](#) (bool b)  
*By default set Sequence & [Item](#) length to Undefined to avoid recomputing length:*
- void [SetRecomputeSequenceLength](#) (bool b)
- void [SetUseVRUN](#) (bool b)  
*When [VR](#)=16bits in explicit but Implicit has a 32bits length, use [VR](#)=UN.*

### Protected Member Functions

- bool [ChangeFMI](#) ()
- bool [ProcessDataSet](#) ([DataSet](#) &ds, [Dicts](#) const &dicts)

### 27.122.1 Detailed Description

[FileExplicitFilter](#) class After changing a file from Implicit to Explicit representation (see [ImageChangeTransferSyntax](#)) one operation is to make sure the [VR](#) of each DICOM attribute are accurate and do match the one from PS 3.6. Indeed when a file is written in Implicit representation, the [VR](#) is not stored directly in the file.

### Warning

changing an implicit dataset to an explicit dataset is NOT a trivial task of simply changing the [VR](#) to the dict one:

- One has to make sure SQ is properly set
- One has to recompute the explicit length SQ
- One has to make sure that [VR](#) is valid for the encoding
- One has to make sure that [VR](#) 16bits can store the original value length

### Examples:

[GenAllVR.cxx](#), and [LargeVRDSExplicit.cxx](#).

## 27.122.2 Constructor & Destructor Documentation

27.122.2.1 `gdcmm::FileExplicitFilter::FileExplicitFilter ( )` `[inline]`

27.122.2.2 `gdcmm::FileExplicitFilter::~~FileExplicitFilter ( )` `[inline]`

## 27.122.3 Member Function Documentation

27.122.3.1 `bool gdcmm::FileExplicitFilter::Change ( )`

Set FMI Transfer Syntax.

Change

### Examples:

[GenAllVR.cxx](#), and [LargeVRDSExplicit.cxx](#).

27.122.3.2 `bool gdcmm::FileExplicitFilter::ChangeFMI ( )` `[protected]`

27.122.3.3 `File& gdcmm::FileExplicitFilter::GetFile ( )` `[inline]`

27.122.3.4 `bool gdcmm::FileExplicitFilter::ProcessDataSet ( DataSet & ds, Dicts const & dicts )` `[protected]`

27.122.3.5 `void gdcmm::FileExplicitFilter::SetChangePrivateTags ( bool b )` `[inline]`

Decide whether or not to [VR](#)ify private tags.

27.122.3.6 `void gdcmm::FileExplicitFilter::SetFile ( const File & f )` `[inline]`

Set/Get [File](#).

### Examples:

[GenAllVR.cxx](#), and [LargeVRDSExplicit.cxx](#).

27.122.3.7 `void gdcmm::FileExplicitFilter::SetRecomputeItemLength ( bool b )`

By default set Sequence & [Item](#) length to Undefined to avoid recomputing length:

27.122.3.8 void gdcm::FileExplicitFilter::SetRecomputeSequenceLength ( bool *b* )

27.122.3.9 void gdcm::FileExplicitFilter::SetUseVRUN ( bool *b* ) [inline]

When [VR](#)=16bits in explicit but Implicit has a 32bits length, use [VR](#)=UN.

The documentation for this class was generated from the following file:


- [gdcmFileExplicitFilter.h](#)

## 27.123 gdcm::FileMetaInformation Class Reference


Class to represent a [File](#) Meta Information.

```
#include <gdcmFileMetaInformation.h>
```

Inheritance diagram for gdcm::FileMetaInformation:

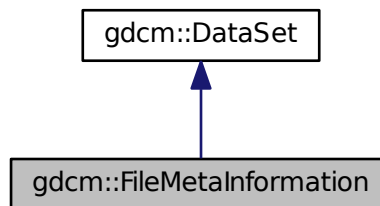

Collaboration diagram for gdcm::FileMetaInformation:

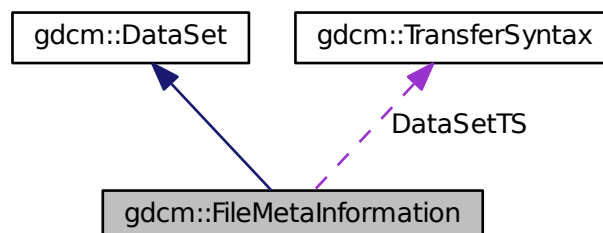

## Public Member Functions

- [FileMetaInformation](#) ()
- [FileMetaInformation](#) ([FileMetaInformation](#) const &fmi)
- void [FillFromDataSet](#) ([DataSet](#) const &ds)
 

*Construct a [FileMetaInformation](#) from an already existing [DataSet](#):*
- const [TransferSyntax](#) & [GetDataSetTransferSyntax](#) () const
- [VL](#) [GetFullLength](#) () const
- [MediaStorage](#) [GetMediaStorage](#) () const
- std::string [GetMediaStorageAsString](#) () const
- [TransferSyntax::NegociatedType](#) [GetMetaInformationTS](#) () const
- const [Preamble](#) & [GetPreamble](#) () const
 

*Get [Preamble](#).*
- [Preamble](#) & [GetPreamble](#) ()
- void [Insert](#) (const [DataElement](#) &de)
- bool [IsValid](#) () const
- std::istream & [Read](#) (std::istream &is)
 

*Read.*
- std::istream & [ReadCompat](#) (std::istream &is)
- void [Replace](#) (const [DataElement](#) &de)
- void [SetDataSetTransferSyntax](#) (const [TransferSyntax](#) &ts)
- void [SetPreamble](#) (const [Preamble](#) &p)
- std::ostream & [Write](#) (std::ostream &os) const
 

*Write.*

## Static Public Member Functions

- static void [AppendImplementationClassUID](#) (const char \*imp)
- static const char \* [GetImplementationClassUID](#) ()
- static const char \* [GetImplementationVersionName](#) ()
- static const char \* [GetSourceApplicationEntityTitle](#) ()
- static void [SetImplementationClassUID](#) (const char \*imp)
 

*Override the GDCM default values:*
- static void [SetImplementationVersionName](#) (const char \*version)
- static void [SetSourceApplicationEntityTitle](#) (const char \*title)

## Protected Member Functions

- void [ComputeDataSetMediaStorageSOPClass](#) ()
- void [ComputeDataSetTransferSyntax](#) ()
- void [Default](#) ()
- template<typename TSwap >  
std::istream & [ReadCompatInternal](#) (std::istream &is)

## Static Protected Member Functions

- static const char \* [GetFileMetaInformationVersion](#) ()
- static const char \* [GetGDCMImplementationClassUID](#) ()
- static const char \* [GetGDCMImplementationVersionName](#) ()
- static const char \* [GetGDCMSourceApplicationEntityTitle](#) ()



## Protected Attributes

- [MediaStorage::MSType DataSetMS](#)
- [TransferSyntax DataSetTS](#)
- [TransferSyntax::NegociatedType MetaInformationTS](#)

## Friends

- `std::ostream & operator<< (std::ostream &_os, const FileMetaInformation &_val)`

## Additional Inherited Members

### 27.123.1 Detailed Description

Class to represent a [File](#) Meta Information.

[FileMetaInformation](#) is a Explicit Structured Set. Whenever the file contains an [ImplicitDataElement DataSet](#), a conversion will take place.

Definition: The [File](#) Meta Information includes identifying information on the encapsulated Data Set. This header consists of a 128 byte [File Preamble](#), followed by a 4 byte DICOM prefix, followed by the [File](#) Meta Elements shown in [Table 7.1-1](#). This header shall be present in every DICOM file.

See also

[Writer Reader](#)

Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [LargeVRDSEExplicit.cxx](#), and [ReadAndDumpDICOMDIR.cxx](#).

### 27.123.2 Constructor & Destructor Documentation

27.123.2.1 `gdcm::FileMetaInformation::FileMetaInformation ( )` `[inline]`

27.123.2.2 `gdcm::FileMetaInformation::FileMetaInformation ( FileMetaInformation const & fmi )` `[inline]`

References [DataSetMS](#), [DataSetTS](#), and [MetaInformationTS](#).

### 27.123.3 Member Function Documentation

27.123.3.1 `static void gdcm::FileMetaInformation::AppendImplementationClassUID ( const char * imp )` `[static]`

27.123.3.2 `void gdcm::FileMetaInformation::ComputeDataSetMediaStorageSOPClass ( )` `[protected]`

27.123.3.3 `void gdcm::FileMetaInformation::ComputeDataSetTransferSyntax ( )` `[protected]`

27.123.3.4 `void gdcm::FileMetaInformation::Default ( )` `[protected]`

27.123.3.5 `void gdcm::FileMetaInformation::FillFromDataSet ( DataSet const & ds )`

Construct a [FileMetaInformation](#) from an already existing [DataSet](#):

27.123.3.6 `const TransferSyntax& gdcm::FileMetaInformation::GetDataSetTransferSyntax ( ) const` `[inline]`

Examples:

[GetJPEGSamplePrecision.cxx](#), and [MergeTwoFiles.cxx](#).

27.123.3.7 `static const char* gdcm::FileMetaInformation::GetFileMetaInformationVersion ( )` `[static]`, `[protected]`

27.123.3.8 `VL gdcm::FileMetaInformation::GetFullLength ( ) const` `[inline]`

References `gdcm::VL::GetLength()`.

27.123.3.9 `static const char* gdcm::FileMetaInformation::GetGDCMImplementationClassUID ( )` `[static]`, `[protected]`

27.123.3.10 `static const char* gdcm::FileMetaInformation::GetGDCMImplementationVersionName ( )` `[static]`, `[protected]`

27.123.3.11 `static const char* gdcm::FileMetaInformation::GetGDCMSourceApplicationEntityTitle ( )` `[static]`, `[protected]`

27.123.3.12 `static const char* gdcm::FileMetaInformation::GetImplementationClassUID ( )` `[static]`

27.123.3.13 `static const char* gdcm::FileMetaInformation::GetImplementationVersionName ( )` `[static]`

27.123.3.14 `MediaStorage gdcm::FileMetaInformation::GetMediaStorage ( ) const`

27.123.3.15 `std::string gdcm::FileMetaInformation::GetMediaStorageAsString ( ) const`

27.123.3.16 `TransferSyntax::NegociatedType gdcm::FileMetaInformation::GetMetaInformationTS ( ) const` `[inline]`

27.123.3.17 `const Preamble& gdcm::FileMetaInformation::GetPreamble ( ) const` `[inline]`

Get [Preamble](#).

Referenced by `gdcm::operator<<()`.

27.123.3.18 `Preamble& gdcm::FileMetaInformation::GetPreamble ( )` `[inline]`

27.123.3.19 `static const char* gdcm::FileMetaInformation::GetSourceApplicationEntityTitle ( )` `[static]`

27.123.3.20 `void gdcm::FileMetaInformation::Insert ( const DataElement & de )` `[inline]`

References `gdcmErrorMacro`, `gdcm::Tag::GetGroup()`, and `gdcm::DataElement::GetTag()`.

27.123.3.21 `bool gdcm::FileMetaInformation::IsValid ( ) const` `[inline]`

27.123.3.22 `std::istream& gdcm::FileMetaInformation::Read ( std::istream & is )`

Read.

27.123.3.23 `std::istream& gdcmm::FileMetaInformation::ReadCompat ( std::istream & is )`

27.123.3.24 `template<typename TSwap > std::istream& gdcmm::FileMetaInformation::ReadCompatInternal ( std::istream & is )`  
[protected]

27.123.3.25 `void gdcmm::FileMetaInformation::Replace ( const DataElement & de )` [inline]

Examples:

[LargeVRDSExplicit.cxx](#).

References `gdcmm::DataElement::GetTag()`.

27.123.3.26 `void gdcmm::FileMetaInformation::SetDataSetTransferSyntax ( const TransferSyntax & ts )`

Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [LargeVRDSExplicit.cxx](#), [MpegVideoInfo.cs](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [rle2img.cxx](#), and [StreamImageReaderTest.cxx](#).

27.123.3.27 `static void gdcmm::FileMetaInformation::SetImplementationClassUID ( const char * imp )` [static]

Override the GDCM default values:

27.123.3.28 `static void gdcmm::FileMetaInformation::SetImplementationVersionName ( const char * version )` [static]

27.123.3.29 `void gdcmm::FileMetaInformation::SetPreamble ( const Preamble & p )` [inline]

27.123.3.30 `static void gdcmm::FileMetaInformation::SetSourceApplicationEntityTitle ( const char * title )` [static]

Examples:

[ClinicalTrialIdentificationWorkflow.cs](#), [FixJAIBugJPEGLS.cxx](#), [GenerateDICOMDIR.cs](#), [ReformatFile.cs](#), and [StandardizeFiles.cs](#).

27.123.3.31 `std::ostream& gdcmm::FileMetaInformation::Write ( std::ostream & os ) const`

Write.

## 27.123.4 Friends And Related Function Documentation

27.123.4.1 `std::ostream& operator<< ( std::ostream & _os, const FileMetaInformation & _val )` [friend]

## 27.123.5 Member Data Documentation

27.123.5.1 `MediaStorage::MSType gdcmm::FileMetaInformation::DataSetMS` [protected]

Referenced by `FileMetaInformation()`.

### 27.123.5.2 TransferSyntax gdcm::FileMetaInformation::DataSetTS [protected]

Referenced by FileMetaInformation().

### 27.123.5.3 TransferSyntax::NegociatedType gdcm::FileMetaInformation::MetaInformationTS [protected]

Referenced by FileMetaInformation().

The documentation for this class was generated from the following file:

- [gdcmFileMetaInformation.h](#)

## 27.124 gdcm::Filename Class Reference

Class to manipulate file name's.

```
#include <gdcmFilename.h>
```

### Public Member Functions

- [Filename](#) (const char \*filename="")
- bool [EndWith](#) (const char ending[]) const  
*Does the filename ends with a particular string ?*
- const char \* [GetExtension](#) ()  
*return only the extension part of a filename*
- const char \* [GetFileName](#) () const  
*Return the full filename.*
- const char \* [GetName](#) ()  
*return only the name part of a filename*
- const char \* [GetPath](#) ()  
*Return only the path component of a filename.*
- bool [IsEmpty](#) () const  
*return whether the filename is empty*
- bool [IsIdentical](#) ([Filename](#) const &fn) const
- [operator const char \\*](#) () const
- const char \* [ToUnixSlashes](#) ()  
*Convert backslash (windows style) to UNIX style slash.*
- const char \* [ToWindowsSlashes](#) ()  
*Convert foward slash (UNIX style) to windows style slash.*

### Static Public Member Functions

- static const char \* [Join](#) (const char \*path, const char \*filename)

### 27.124.1 Detailed Description

Class to manipulate file name's.

#### Note

OS independant representation of a filename (to query path, name and extension from a filename)

### 27.124.2 Constructor & Destructor Documentation

27.124.2.1 `gdcm::Filename::Filename ( const char * filename = " " ) [inline]`

### 27.124.3 Member Function Documentation

27.124.3.1 `bool gdcm::Filename::EndWith ( const char ending[ ] ) const`

Does the filename ends with a particular string ?

27.124.3.2 `const char* gdcm::Filename::GetExtension ( )`

return only the extension part of a filename

27.124.3.3 `const char* gdcm::Filename::GetFileName ( ) const [inline]`

Return the full filename.

27.124.3.4 `const char* gdcm::Filename::GetName ( )`

return only the name part of a filename

27.124.3.5 `const char* gdcm::Filename::GetPath ( )`

Return only the path component of a filename.

#### Examples:

[ClinicalTrialIdentificationWorkflow.cs](#).

27.124.3.6 `bool gdcm::Filename::IsEmpty ( ) const [inline]`

return whether the filename is empty

27.124.3.7 `bool gdcm::Filename::IsIdentical ( Filename const & fn ) const`

27.124.3.8 `static const char* gdcm::Filename::Join ( const char * path, const char * filename ) [static]`

Join two paths NOT THREAD SAFE

Examples:

[BasicAnonymizer.cs](#), and [ClinicalTrialIdentificationWorkflow.cs](#).

27.124.3.9 `gdcmm::Filename::operator const char * ( ) const` `[inline]`

Simple operator to allow [Filename](#) `myfilename( "...")`; `const char * s = myfilename`;

27.124.3.10 `const char* gdcmm::Filename::ToUnixSlashes ( )`

Convert backslash (windows style) to UNIX style slash.

27.124.3.11 `const char* gdcmm::Filename::ToWindowsSlashes ( )`

Convert forward slash (UNIX style) to windows style slash.

The documentation for this class was generated from the following file:


- [gdcmmFilename.h](#)

## 27.125 `gdcmm::FileNameEvent` Class Reference


[FileNameEvent](#) Special type of event triggered during processing of [FileSet](#).

```
#include <gdcmmFileNameEvent.h>
```

Inheritance diagram for `gdcmm::FileNameEvent`:

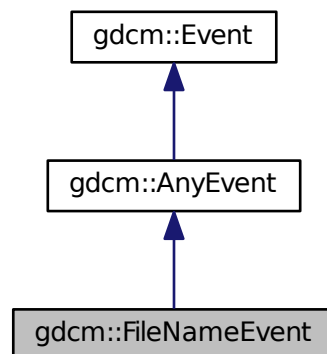

Collaboration diagram for gdcm::FileNameEvent:

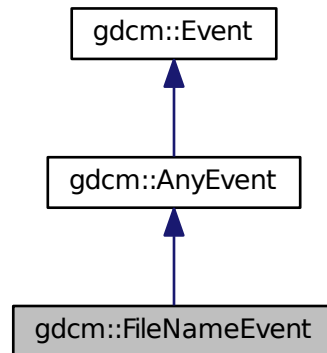

## Public Types

- typedef [FileNameEvent](#) Self
- typedef [AnyEvent](#) Superclass

## Public Member Functions

- [FileNameEvent](#) (const char \*s="")
- [FileNameEvent](#) (const Self &s)
- virtual [~FileNameEvent](#) ()
- virtual bool [CheckEvent](#) (const ::gdcm::Event \*e) const
- virtual const char \* [GetEventName](#) () const
- const char \* [GetFileName](#) () const
- virtual ::gdcm::Event \* [MakeObject](#) () const
- void [SetFileName](#) (const char \*f)

### 27.125.1 Detailed Description

[FileNameEvent](#) Special type of event triggered during processing of [FileSet](#).

See also

[AnyEvent](#)

Examples:

[SimpleScanner.cxx](#).

## 27.125.2 Member Typedef Documentation

27.125.2.1 `typedef FileNameEvent gdcm::FileNameEvent::Self`

27.125.2.2 `typedef AnyEvent gdcm::FileNameEvent::Superclass`

## 27.125.3 Constructor & Destructor Documentation

27.125.3.1 `gdcm::FileNameEvent::FileNameEvent ( const char * s = " " ) [inline]`

27.125.3.2 `virtual gdcm::FileNameEvent::~~FileNameEvent ( ) [inline],[virtual]`

27.125.3.3 `gdcm::FileNameEvent::FileNameEvent ( const Self & s ) [inline]`

## 27.125.4 Member Function Documentation

27.125.4.1 `virtual bool gdcm::FileNameEvent::CheckEvent ( const ::gdcm::Event * e ) const [inline],[virtual]`

27.125.4.2 `virtual const char* gdcm::FileNameEvent::GetEventName ( ) const [inline],[virtual]`

Return the StringName associated with the event.

Implements [gdcm::Event](#).

27.125.4.3 `const char* gdcm::FileNameEvent::GetFileName ( ) const [inline]`

Examples:

[SimpleScanner.cxx](#).

27.125.4.4 `virtual ::gdcm::Event* gdcm::FileNameEvent::MakeObject ( ) const [inline],[virtual]`

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implements [gdcm::Event](#).

27.125.4.5 `void gdcm::FileNameEvent::SetFileName ( const char * f ) [inline]`

The documentation for this class was generated from the following file:

- [gdcmFileNameEvent.h](#)

## 27.126 gdcm::FilenameGenerator Class Reference

[FilenameGenerator](#).

```
#include <gdcmFilenameGenerator.h>
```



## Public Types

- typedef std::vector< [FilenameType](#) > [FileNamesType](#)
- typedef std::string [FilenameType](#)
- typedef FileNamesType::size\_type [SizeType](#)

## Public Member Functions

- [FilenameGenerator](#) ()
- [~FilenameGenerator](#) ()
- bool [Generate](#) ()  
*Generate (return success)*
- const char \* [GetFilename](#) ([SizeType](#) n) const  
*Get a particular filename (call after Generate)*
- [FileNamesType](#) const & [GetFileNames](#) () const  
*Return all filenames.*
- [SizeType](#) [GetNumberOfFileNames](#) () const
- const char \* [GetPattern](#) () const
- const char \* [GetPrefix](#) () const
- void [SetNumberOfFileNames](#) ([SizeType](#) nfiles)  
*Set/Get the number of filenames to generate.*
- void [SetPattern](#) (const char \*pattern)  
*Set/Get pattern.*
- void [SetPrefix](#) (const char \*prefix)  
*Set/Get prefix.*

### 27.126.1 Detailed Description

[FilenameGenerator](#).

class to generate filenames based on a pattern (C-style)

Output will be:

for i = 0, number of filenames: outfilename[i] = prefix + (pattern % i)

where pattern % i means C-style sprintf of Pattern using value 'i'

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

### 27.126.2 Member Typedef Documentation

27.126.2.1 typedef std::vector<[FilenameType](#)> gdcm::FilenameGenerator::FileNamesType

27.126.2.2 typedef std::string gdcm::FilenameGenerator::FilenameType

27.126.2.3 typedef FileNamesType::size\_type gdcm::FilenameGenerator::SizeType

### 27.126.3 Constructor & Destructor Documentation

27.126.3.1 `gdcm::FilenameGenerator::FilenameGenerator ( )` `[inline]`

27.126.3.2 `gdcm::FilenameGenerator::~~FilenameGenerator ( )` `[inline]`

## 27.126.4 Member Function Documentation

27.126.4.1 `bool gdcm::FilenameGenerator::Generate ( )`

Generate (return success)

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.126.4.2 `const char* gdcm::FilenameGenerator::GetFilename ( SizeType n ) const`

Get a particular filename (call after Generate)

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.126.4.3 `FilenameType const& gdcm::FilenameGenerator::GetFilenames ( ) const` `[inline]`

Return all filenames.

27.126.4.4 `SizeType gdcm::FilenameGenerator::GetNumberOfFilenames ( ) const`

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.126.4.5 `const char* gdcm::FilenameGenerator::GetPattern ( ) const` `[inline]`

27.126.4.6 `const char* gdcm::FilenameGenerator::GetPrefix ( ) const` `[inline]`

27.126.4.7 `void gdcm::FilenameGenerator::SetNumberOfFilenames ( SizeType nfiles )`

Set/Get the number of filenames to generate.

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.126.4.8 `void gdcm::FilenameGenerator::SetPattern ( const char * pattern )` `[inline]`

Set/Get pattern.

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.126.4.9 void gdcm::FilenameGenerator::SetPrefix ( const char \* *prefix* ) [inline]

Set/Get prefix.

The documentation for this class was generated from the following file:

- [gdcmFilenameGenerator.h](#)

## 27.127 gdcm::FileSet Class Reference

File-set: A File-set is a collection of DICOM Files (and possibly non-DICOM Files) that share a common naming space within which [File](#) IDs are unique.

```
#include <gdcmFileSet.h>
```

### Public Types

- typedef std::vector< [FileType](#) > [FilesType](#)
- typedef std::string [FileType](#)

### Public Member Functions

- [FileSet](#) ()
- void [AddFile](#) ([File](#) const &)
- bool [AddFile](#) (const char \*filename)
- [FilesType](#) const & [GetFiles](#) () const
- void [SetFiles](#) ([FilesType](#) const &files)

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [FileSet](#) &d)

### 27.127.1 Detailed Description

File-set: A File-set is a collection of DICOM Files (and possibly non-DICOM Files) that share a common naming space within which [File](#) IDs are unique.

### 27.127.2 Member Typedef Documentation

27.127.2.1 typedef std::vector<[FileType](#)> gdcm::FileSet::FilesType

27.127.2.2 typedef std::string gdcm::FileSet::FileType

### 27.127.3 Constructor & Destructor Documentation

27.127.3.1 gdcm::FileSet::FileSet ( ) [inline]

### 27.127.4 Member Function Documentation

27.127.4.1 `void gdcM::FileSet::AddFile ( File const & ) [inline]`

**Deprecated** . Does nothing

27.127.4.2 `bool gdcM::FileSet::AddFile ( const char * filename )`

Add a file 'filename' to the list of files. Return true on success, false in case filename could not be found on system.

27.127.4.3 `FileType const& gdcM::FileSet::GetFiles ( ) const [inline]`

27.127.4.4 `void gdcM::FileSet::SetFiles ( FileType const & files )`

## 27.127.5 Friends And Related Function Documentation

27.127.5.1 `std::ostream& operator<< ( std::ostream & _os, const FileSet & d ) [friend]`

The documentation for this class was generated from the following file:


- [gdcMFileSet.h](#)

## 27.128 gdcM::FileStreamer Class Reference


**FileStreamer** This class let a user create a massive DICOM [DataSet](#) from a template DICOM file, by appending chunks of data.

```
#include <gdcMFileStreamer.h>
```

Inheritance diagram for gdcM::FileStreamer:

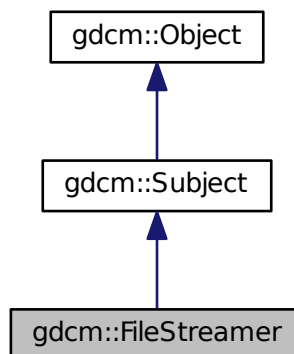

Collaboration diagram for gdcm::FileStreamer:

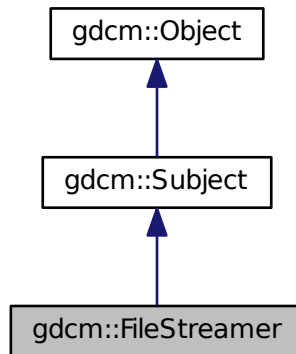

## Public Member Functions

- [FileStreamer](#) ()
- [~FileStreamer](#) ()
- bool [AppendToDataElement](#) (const [Tag](#) &t, const char \*array, size\_t len)  
*Append to previously started Tag t.*
- bool [AppendToGroupDataElement](#) (const [PrivateTag](#) &pt, const char \*array, size\_t len)  
*Append to previously started private creator.*
- bool [CheckDataElement](#) (const [Tag](#) &t)
- void [CheckTemplateFileName](#) (bool check)
- bool [ReserveDataElement](#) (size\_t len)
- bool [ReserveGroupDataElement](#) (unsigned short ndataelement)
- void [SetOutputFileName](#) (const char \*filename\_native)  
*Set output filename (target file)*
- void [SetTemplateFileName](#) (const char \*filename\_native)  
*Set input DICOM template filename.*
- bool [StartDataElement](#) (const [Tag](#) &t)
- bool [StartGroupDataElement](#) (const [PrivateTag](#) &pt, size\_t maxsize=0, uint8\_t startoffset=0)
- bool [StopDataElement](#) (const [Tag](#) &t)  
*Stop appending to tag t. This will compute the proper attribute length.*
- bool [StopGroupDataElement](#) (const [PrivateTag](#) &pt)  
*Stop appending to private creator.*

## Static Public Member Functions

- static [SmartPointer](#)< [FileStreamer](#) > [New](#) ()  
*for wrapped language: instantiate a reference counted object*

## Additional Inherited Members

### 27.128.1 Detailed Description

**FileStreamer** This class let a user create a massive DICOM [DataSet](#) from a template DICOM file, by appending chunks of data.

This class support two mode of operation:

1. Creating a single [DataElement](#) by appending chunk after chunk of data.
2. Creating a set of [DataElement](#) within the same group, using a private creator for start. New [DataElement](#) are added any time the user defined maximum size for data element is reached.

#### Warning

any existing [DataElement](#) is removed, pick carefully which [DataElement](#) to add.

### 27.128.2 Constructor & Destructor Documentation

27.128.2.1 `gdcm::FileStreamer::FileStreamer ( )`

27.128.2.2 `gdcm::FileStreamer::~~FileStreamer ( )`

### 27.128.3 Member Function Documentation

27.128.3.1 `bool gdcm::FileStreamer::AppendToDataElement ( const Tag & t, const char * array, size_t len )`

Append to previously started [Tag](#) t.

27.128.3.2 `bool gdcm::FileStreamer::AppendToGroupDataElement ( const PrivateTag & pt, const char * array, size_t len )`

Append to previously started private creator.

27.128.3.3 `bool gdcm::FileStreamer::CheckDataElement ( const Tag & t )`

Decide to check the Data [Element](#) to be written (default: off) The implementation has default strategy for checking validity of [DataElement](#). Currently it only support checking for the following tags:

- (7fe0,0010) Pixel Data

27.128.3.4 `void gdcm::FileStreamer::CheckTemplateFileName ( bool check )`

Instead of simply blindly copying the input DICOM Template file, GDCM will be used to check the input file, and correct any issues recognized within the file. Only use if you do not have control over the input template file.

27.128.3.5 `static SmartPointer<FileStreamer> gdcm::FileStreamer::New ( ) [inline], [static]`

for wrapped language: instantiate a reference counted object

27.128.3.6 `bool gdcm::FileStreamer::ReserveDataElement ( size_t len )`

Add a hint on the final size of the dataelement. When optimally chosen, this reduce the number of file in-place copying. Should be called before StartDataElement

27.128.3.7 `bool gdcm::FileStreamer::ReserveGroupDataElement ( unsigned short ndataelement )`

Optimisation: pre-allocate the number of dataelement within the private group (ndataelement <= 256). Should be called before StartGroupDataElement

27.128.3.8 `void gdcm::FileStreamer::SetOutputFileName ( const char * filename_native )`

Set output filename (target file)

27.128.3.9 `void gdcm::FileStreamer::SetTemplateFileName ( const char * filename_native )`

Set input DICOM template filename.

Examples:

[FileStreaming.cs](#).

27.128.3.10 `bool gdcm::FileStreamer::StartDataElement ( const Tag & t )`

Start Single Data [Element](#) Operation This will delete any existing [Tag](#) t. Need to call it only once.

27.128.3.11 `bool gdcm::FileStreamer::StartGroupDataElement ( const PrivateTag & pt, size_t maxsize = 0, uint8_t startoffset = 0 )`

Start Private Group (multiple [DataElement](#)) Operation. Each newly added [DataElement](#) will have a length lower than

Parameters

<i>maxsize</i>	. When not specified, maxsize is set to maximum size allowed by DICOM (= 2 <sup>32</sup> ). startoffset can be used to specify the very first element you want to start with (instead of the first possible). <a href="#">Value</a> should be in [0x0, 0xff] This will find the first available private creator.
----------------	--

27.128.3.12 `bool gdcm::FileStreamer::StopDataElement ( const Tag & t )`

Stop appending to tag t. This will compute the proper attribute length.

27.128.3.13 `bool gdcm::FileStreamer::StopGroupDataElement ( const PrivateTag & pt )`

Stop appending to private creator.

The documentation for this class was generated from the following file:


- [gdcmFileStreamer.h](#)

## 27.129 gdcm::FileWithName Class Reference


[FileWithName.](#)

```
#include <gdcmSerieHelper.h>
```

Inheritance diagram for `gdcm::FileWithName`:

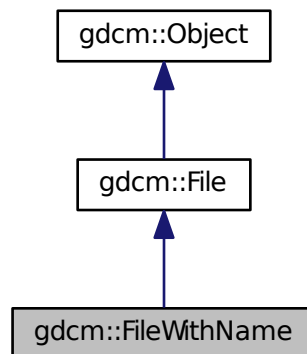

Collaboration diagram for `gdcm::FileWithName`:

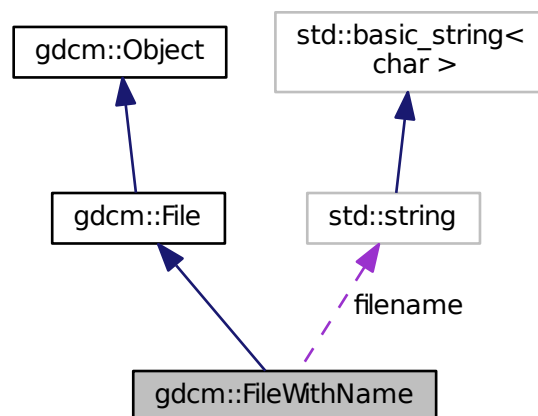

### Public Member Functions

- [FileWithName](#) (`File` &`f`)



## Public Attributes

- `std::string` [filename](#)

## Additional Inherited Members

### 27.129.1 Detailed Description

[FileWithName](#).

Backward only class do not use in newer code

### 27.129.2 Constructor & Destructor Documentation

27.129.2.1 `gdcm::FileWithName::FileWithName ( File & f )` `[inline]`

### 27.129.3 Member Data Documentation

27.129.3.1 `std::string` `gdcm::FileWithName::filename`

The documentation for this class was generated from the following file:


- [gdcmSerieHelper.h](#)

## 27.130 gdcm::FindPatientRootQuery Class Reference


PatientRootQuery contains: the class which will produce a dataset for c-find with patient root.

```
#include <gdcmFindPatientRootQuery.h>
```

Inheritance diagram for `gdcm::FindPatientRootQuery`:

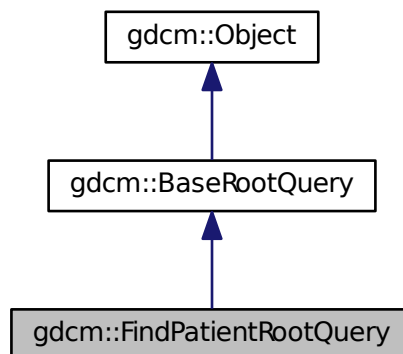

Collaboration diagram for `gdcm::FindPatientRootQuery`:

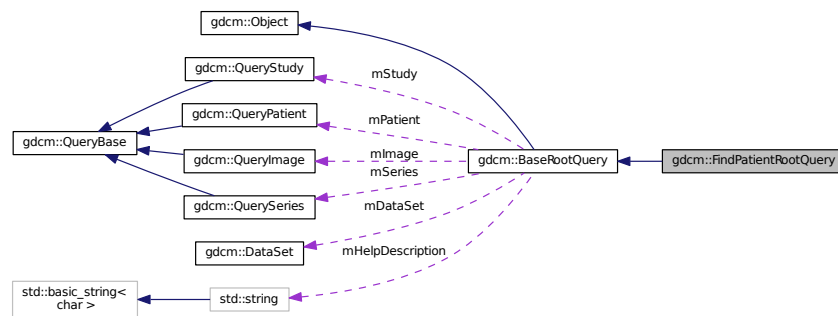

## Public Member Functions

- [FindPatientRootQuery](#) ()
- [UIDs::TSName GetAbstractSyntaxUID](#) () const
- `std::vector< Tag > GetTagListByLevel` (const [EQueryLevel](#) &inQueryLevel)
- void [InitializeDataSet](#) (const [EQueryLevel](#) &inQueryLevel)
- bool [ValidateQuery](#) (bool inStrict=true) const

## Friends

- class [QueryFactory](#)

## Additional Inherited Members

### 27.130.1 Detailed Description

PatientRootQuery contains: the class which will produce a dataset for c-find with patient root.

### 27.130.2 Constructor & Destructor Documentation

27.130.2.1 `gdcm::FindPatientRootQuery::FindPatientRootQuery ( )`

### 27.130.3 Member Function Documentation

27.130.3.1 `UIDs::TSName gdcm::FindPatientRootQuery::GetAbstractSyntaxUID ( ) const` `[virtual]`

Implements [gdcm::BaseRootQuery](#).

27.130.3.2 `std::vector<Tag> gdcm::FindPatientRootQuery::GetTagListByLevel ( const EQueryLevel & inQueryLevel )` `[virtual]`

this function will return all tags at a given query level, so that they maybe selected for searching. The boolean forFind is true if the query is a find query, or false for a move query.

Implements [gdcm::BaseRootQuery](#).

27.130.3.3 void gdcm::FindPatientRootQuery::InitializeDataSet ( const EQueryLevel & inQueryLevel ) [virtual]

this function sets tag 8,52 to the appropriate value based on query level also fills in the right unique tags, as per the standard's requirements should allow for connection with dcmTk

Implements [gdcm::BaseRootQuery](#).

27.130.3.4 bool gdcm::FindPatientRootQuery::ValidateQuery ( bool inStrict = true ) const [virtual]

have to be able to ensure that 0x8,0x52 is set (which will be true if InitializeDataSet is called...) that the level is appropriate (ie, not setting PATIENT for a study query that the tags in the query match the right level (either required, unique, optional) by default, this function checks to see if the query is for finding, which is more permissive than for moving. For moving, only the unique tags are allowed. 10 Jan 2011: adding in the 'strict' mode. according to the standard (at least, how I've read it), only tags for a particular level should be allowed in a particular query (ie, just series level tags in a series level query). However, it seems that dcm4chee doesn't share that interpretation. So, if 'inStrict' is false, then tags from the current level and all higher levels are now considered valid. So, if you're doing a non-strict series-level query, tags from the patient and study level can be passed along as well.

Implements [gdcm::BaseRootQuery](#).

## 27.130.4 Friends And Related Function Documentation

27.130.4.1 friend class QueryFactory [friend]

The documentation for this class was generated from the following file:


- [gdcmFindPatientRootQuery.h](#)

## 27.131 gdcm::FindStudyRootQuery Class Reference


[FindStudyRootQuery](#) contains: the class which will produce a dataset for C-FIND with study root.

```
#include <gdcmFindStudyRootQuery.h>
```

Inheritance diagram for `gdcM::FindStudyRootQuery`:

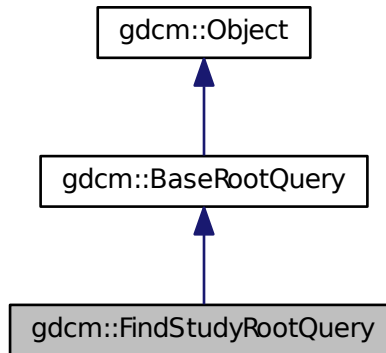

Collaboration diagram for `gdcM::FindStudyRootQuery`:

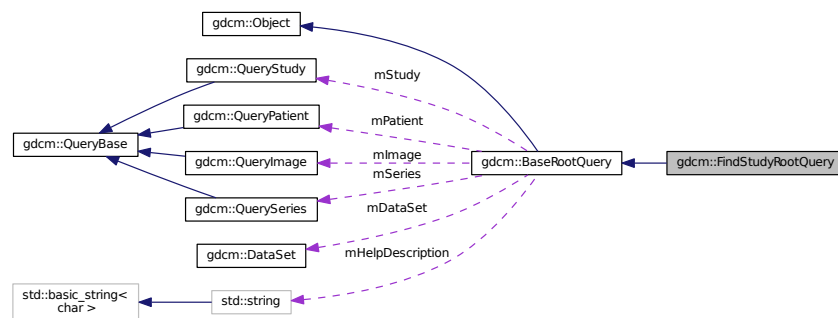

## Public Member Functions

- [FindStudyRootQuery](#) ()
- [UIDs::TSName GetAbstractSyntaxUID](#) () const
- `std::vector< Tag > GetTagListByLevel` (const [EQueryLevel](#) &inQueryLevel)
- void [InitializeDataSet](#) (const [EQueryLevel](#) &inQueryLevel)
- bool [ValidateQuery](#) (bool inStrict=true) const

## Friends

- class [QueryFactory](#)

## Additional Inherited Members

### 27.131.1 Detailed Description

[FindStudyRootQuery](#) contains: the class which will produce a dataset for C-FIND with study root.

### 27.131.2 Constructor & Destructor Documentation

27.131.2.1 `gdcm::FindStudyRootQuery::FindStudyRootQuery ( )`

### 27.131.3 Member Function Documentation

27.131.3.1 `UIDs::TSName gdcm::FindStudyRootQuery::GetAbstractSyntaxUID ( ) const [virtual]`

Implements [gdcm::BaseRootQuery](#).

27.131.3.2 `std::vector<Tag> gdcm::FindStudyRootQuery::GetTagListByLevel ( const EQueryLevel & inQueryLevel ) [virtual]`

this function will return all tags at a given query level, so that they maybe selected for searching. The boolean forFind is true if the query is a find query, or false for a move query.

Implements [gdcm::BaseRootQuery](#).

27.131.3.3 `void gdcm::FindStudyRootQuery::InitializeDataSet ( const EQueryLevel & inQueryLevel ) [virtual]`

this function sets tag 8,52 to the appropriate value based on query level also fills in the right unique tags, as per the standard's requirements should allow for connection with dcmtk

Implements [gdcm::BaseRootQuery](#).

27.131.3.4 `bool gdcm::FindStudyRootQuery::ValidateQuery ( bool inStrict=true ) const [virtual]`

have to be able to ensure that (0008,0052) is set that the level is appropriate (ie, not setting PATIENT for a study query that the tags in the query match the right level (either required, unique, optional)

Implements [gdcm::BaseRootQuery](#).

### 27.131.4 Friends And Related Function Documentation

27.131.4.1 `friend class QueryFactory [friend]`

The documentation for this class was generated from the following file:


- [gdcmFindStudyRootQuery.h](#)

## 27.132 gdcm::Fragment Class Reference


Class to represent a [Fragment](#).

```
#include <gdcmFragment.h>
```

Inheritance diagram for `gdcm::Fragment`:

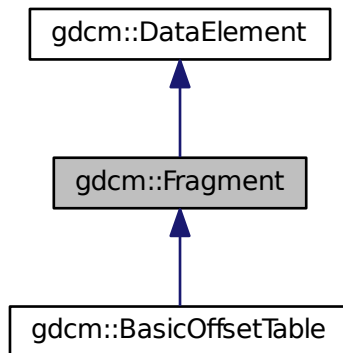

Collaboration diagram for `gdcm::Fragment`:

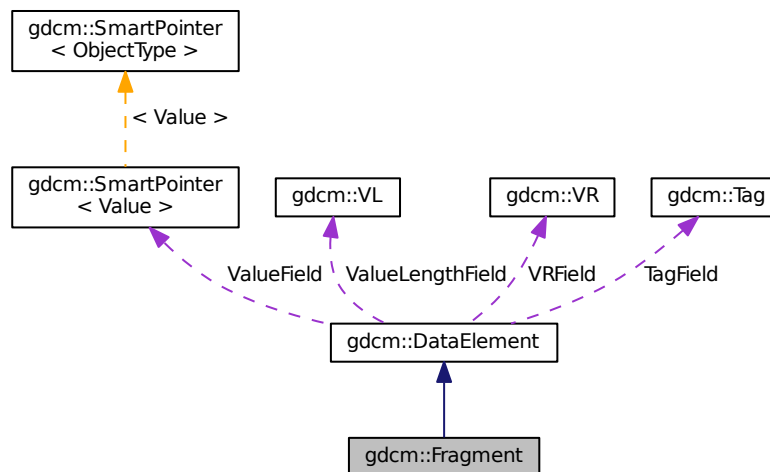

## Public Member Functions

- [Fragment](#) ()
- [VL ComputeLength](#) () const
- [VL GetLength](#) () const

- `template<typename TSwap >`  
`std::istream & Read (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & ReadBacktrack (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & ReadPreValue (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & ReadValue (std::istream &is)`
- `template<typename TSwap >`  
`std::ostream & Write (std::ostream &os) const`

## Friends

- `std::ostream & operator<< (std::ostream &os, const Fragment &val)`

## Additional Inherited Members

### 27.132.1 Detailed Description

Class to represent a [Fragment](#).

Examples:

[FixBrokenJ2K.cxx](#), and [FixJAIBugJPEGLS.cxx](#).

### 27.132.2 Constructor & Destructor Documentation

27.132.2.1 `gdcm::Fragment::Fragment ( )` `[inline]`

### 27.132.3 Member Function Documentation

27.132.3.1 `VL gdcm::Fragment::ComputeLength ( )` `const`

27.132.3.2 `VL gdcm::Fragment::GetLength ( )` `const`

27.132.3.3 `template<typename TSwap > std::istream& gdcm::Fragment::Read ( std::istream & is )` `[inline]`

Referenced by `gdcm::SequenceOfFragments::ReadValue()`.

27.132.3.4 `template<typename TSwap > std::istream& gdcm::Fragment::ReadBacktrack ( std::istream & is )` `[inline]`

References `gdcmErrorMacro`, `gdcmWarningMacro`, and `gdcm::ParseException::SetLastElement()`.

Referenced by `gdcm::SequenceOfFragments::ReadValue()`.

27.132.3.5 `template<typename TSwap > std::istream& gdcm::Fragment::ReadPreValue ( std::istream & is )` `[inline]`

27.132.3.6 `template<typename TSwap > std::istream& gdcm::Fragment::ReadValue ( std::istream & is )` `[inline]`

References `gdcmWarningMacro`, and `gdcm::ParseException::SetLastElement()`.

27.132.3.7 `template<typename TSwap > std::ostream& gdcmm::Fragment::Write ( std::ostream & os ) const` `[inline]`

References `gdcmm::ByteValue::ComputeLength()`, `gdcmm::ByteValue::GetLength()`, `gdcmm::VL::Write()`, and `gdcmm::ByteValue::Write()`.

## 27.132.4 Friends And Related Function Documentation

27.132.4.1 `std::ostream& operator<< ( std::ostream & os, const Fragment & val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmmFragment.h](#)

## 27.133 gdcmm::Global Class Reference

[Global](#).

```
#include <gdcmmGlobal.h>
```

### Public Member Functions

- [Global](#) ()
- [~Global](#) ()
- bool [Append](#) (const char \*path)
- [Defs](#) const & [GetDefs](#) () const
- [Dicts](#) const & [GetDicts](#) () const
- [Dicts](#) & [GetDicts](#) ()
- bool [LoadResourcesFiles](#) ()
- bool [Prepend](#) (const char \*path)

### Static Public Member Functions

- static [Global](#) & [GetInstance](#) ()  
*return the singleton instance*

### Protected Member Functions

- const char \* [Locate](#) (const char \*resfile) const  
*Locate a resource file.*

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Global](#) &g)



### 27.133.1 Detailed Description

[Global](#).

Note

[Global](#) should be included in any translation unit that will use [Dict](#) or that implements the singleton pattern. It makes sure that the [Dict](#) singleton is created before and destroyed after all other singletons in GDCM.

Examples:

[GenAllVR.cxx](#), [GenerateStandardSOPClasses.cxx](#), [GenFakeIdentifyFile.cxx](#), [PublicDict.cxx](#), [ReadAndPrintAttributes.cxx](#), and [TraverseModules.cxx](#).

### 27.133.2 Constructor & Destructor Documentation

27.133.2.1 `gdcmm::Global::Global ( )`

27.133.2.2 `gdcmm::Global::~~Global ( )`

### 27.133.3 Member Function Documentation

27.133.3.1 `bool gdcmm::Global::Append ( const char * path )`

Append path at the end of the path list

Warning

not thread safe !

27.133.3.2 `Defs const& gdcmm::Global::GetDefs ( ) const`

retrieve the default/internal (Part 3) You need to explicitly call LoadResourcesFiles before

Examples:

[GenerateStandardSOPClasses.cxx](#), and [TraverseModules.cxx](#).

27.133.3.3 `Dicts const& gdcmm::Global::GetDicts ( ) const`

retrieve the default/internal dicts (Part 6) This dict is filled up at load time

Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [MrProtocol.cxx](#), [PublicDict.cxx](#), [ReadAndPrintAttributes.cxx](#), and [TraverseModules.cxx](#).

27.133.3.4 `Dicts& gdcmm::Global::GetDicts ( )`

27.133.3.5 `static Global& gdcmm::Global::GetInstance ( ) [static]`

return the singleton instance

**Examples:**

[BasicAnonymizer.cs](#), [ClinicalTrialIdentificationWorkflow.cs](#), [GenAllVR.cxx](#), [GenerateStandardSOPClasses.cxx](#), [GenFakeIdentifyFile.cxx](#), [MrProtocol.cxx](#), [PublicDict.cxx](#), [ReadAndPrintAttributes.cxx](#), and [TraverseModules.cxx](#).

**27.133.3.6 bool gdcm::Global::LoadResourcesFiles ( )**

Load all internal XML files, resource path need to have been set before calling this member function (see [Append/Prepend members func](#))

**Warning**

not thread safe !

**Examples:**

[GenerateStandardSOPClasses.cxx](#), and [TraverseModules.cxx](#).

**27.133.3.7 const char\* gdcm::Global::Locate ( const char \* *resfile* ) const** [protected]

Locate a resource file.

**27.133.3.8 bool gdcm::Global::Prepend ( const char \* *path* )**

Prepend path at the beginning of the path list

**Warning**

not thread safe !

**27.133.4 Friends And Related Function Documentation****27.133.4.1 std::ostream& operator<< ( std::ostream & *\_os*, const Global & *g* )** [friend]

The documentation for this class was generated from the following file:

- [gdcmGlobal.h](#)

**27.134 gdcm::GroupDict Class Reference**

Class to represent the mapping from group number to its abbreviation and name.

```
#include <gdcmGroupDict.h>
```

**Public Types**

- typedef std::vector< std::string > [GroupStringVector](#)

## Public Member Functions

- [GroupDict](#) ()
- [~GroupDict](#) ()
- std::string const & [GetAbbreviation](#) (uint16\_t num) const
- std::string const & [GetName](#) (uint16\_t num) const
- size\_t [Size](#) () const

## Protected Member Functions

- void [Add](#) (std::string const &abbreviation, std::string const &name)
- void [Insert](#) (uint16\_t num, std::string const &abbreviation, std::string const &name)

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [GroupDict](#) &\_val)

### 27.134.1 Detailed Description

Class to represent the mapping from group number to its abbreviation and name.

#### Note

Should I rewrite this class to use a std::map instead of std::vector for problem of memory consumption ?

### 27.134.2 Member Typedef Documentation

27.134.2.1 `typedef std::vector<std::string> gdcm::GroupDict::GroupStringVector`

### 27.134.3 Constructor & Destructor Documentation

27.134.3.1 `gdcm::GroupDict::GroupDict ( )` `[inline]`

27.134.3.2 `gdcm::GroupDict::~~GroupDict ( )` `[inline]`

### 27.134.4 Member Function Documentation

27.134.4.1 `void gdcm::GroupDict::Add ( std::string const & abbreviation, std::string const & name )` `[protected]`

27.134.4.2 `std::string const& gdcm::GroupDict::GetAbbreviation ( uint16_t num ) const`

Referenced by `gdcm::operator<<()`.

27.134.4.3 `std::string const& gdcm::GroupDict::GetName ( uint16_t num ) const`

Referenced by `gdcm::operator<<()`.

27.134.4.4 void gdcM::GroupDict::Insert ( uint16\_t num, std::string const & abbreviation, std::string const & name )  
[protected]

27.134.4.5 size\_t gdcM::GroupDict::Size ( ) const [inline]

Referenced by gdcM::operator<<().

## 27.134.5 Friends And Related Function Documentation

27.134.5.1 std::ostream& operator<< ( std::ostream & \_os, const GroupDict & \_val ) [friend]

The documentation for this class was generated from the following file:

- [gdcMGroupDict.h](#)

## 27.135 gdcM::IconImageFilter Class Reference

**IconImageFilter** This filter will extract icons from a [File](#) This filter will loop over all known sequence (public and private) that may contains an IconImage and retrieve them. The filter will fails with a value of false if no icon can be found Since it handle both public and private icon type, one should not assume the icon is in uncompress form, some private vendor store private icon in JPEG8/JPEG12.

```
#include <gdcMIconImageFilter.h>
```

### Public Member Functions

- [IconImageFilter](#) ()
- [~IconImageFilter](#) ()
- bool [Extract](#) ()  
*Extract all Icon found in [File](#).*
- [File](#) & [GetFile](#) ()
- const [File](#) & [GetFile](#) () const
- [IconImage](#) & [GetIconImage](#) (unsigned int i) const
- unsigned int [GetNumberOfIconImages](#) () const  
*Retrieve extract IconImage (need to call Extract first)*
- void [SetFile](#) (const [File](#) &f)  
*Set/Get [File](#).*

### Protected Member Functions

- void [ExtractIconImages](#) ()
- void [ExtractVeprolIconImages](#) ()

### 27.135.1 Detailed Description

**IconImageFilter** This filter will extract icons from a [File](#) This filter will loop over all known sequence (public and private) that may contains an IconImage and retrieve them. The filter will fails with a value of false if no icon can be found Since

it handle both public and private icon type, one should not assume the icon is in uncompress form, some private vendor store private icon in JPEG8/JPEG12.

Implementation details: This filter supports the following Icons:

- (0088,0200) Icon [Image](#) Sequence
- (0009,10,GEIIS) GE IIS Thumbnail Sequence
- (6003,10,GEMS\_Ultrasound\_ImageGroup\_001) GEMS [Image](#) Thumbnail Sequence
- (0055,30,VEPRO VIF 3.0 DATA) Icon Data
- (0055,30,VEPRO VIM 5.0 DATA) ICONDATA2

#### Warning

the icon stored in those private attribute do not conform to definition of Icon [Image](#) Sequence (do not simply copy/paste). For example some private icon can be expressed as 12bits pixel, while the DICOM standard only allow 8bits icons.

#### See also

[ImageReader](#)

#### Examples:

[ExtractIconFromFile.cxx](#).

### 27.135.2 Constructor & Destructor Documentation

27.135.2.1 `gdcm::IconImageFilter::IconImageFilter ( )`

27.135.2.2 `gdcm::IconImageFilter::~~IconImageFilter ( )`

### 27.135.3 Member Function Documentation

27.135.3.1 `bool gdcm::IconImageFilter::Extract ( )`

Extract all Icon found in [File](#).

#### Examples:

[ExtractIconFromFile.cxx](#).

27.135.3.2 `void gdcm::IconImageFilter::ExtractIconImages ( )` `[protected]`

27.135.3.3 `void gdcm::IconImageFilter::ExtractVeprolconImages ( )` `[protected]`

27.135.3.4 `File& gdcm::IconImageFilter::GetFile ( )` `[inline]`

27.135.3.5 `const File& gdcm::IconImageFilter::GetFile ( ) const` `[inline]`

27.135.3.6 `IconImage& gdcm::IconImageFilter::GetIconImage ( unsigned int i ) const`

Examples:

[ExtractIconFromFile.cxx](#).

27.135.3.7 `unsigned int gdcm::IconImageFilter::GetNumberOfIconImages ( ) const`

Retrieve extract IconImage (need to call Extract first)

Examples:

[ExtractIconFromFile.cxx](#).

27.135.3.8 `void gdcm::IconImageFilter::SetFile ( const File & f ) [inline]`

Set/Get [File](#).

Examples:

[ExtractIconFromFile.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmIconImageFilter.h](#)

## 27.136 gdcm::IconImageGenerator Class Reference

[IconImageGenerator](#) This filter will generate a valid Icon from the Pixel Data element (an instance of [Pixmap](#)). To generate a valid Icon, one is only allowed the following Photometric Interpretation:

```
#include <gdcmIconImageGenerator.h>
```

### Public Member Functions

- [IconImageGenerator](#) ()
- [~IconImageGenerator](#) ()
- void [AutoPixelMinMax](#) (bool b)
- void [ConvertRGBToPaletteColor](#) (bool b)
- bool [Generate](#) ()  
*Generate Icon.*
- const [IconImage](#) & [GetIconImage](#) () const  
*Retrieve generated Icon.*
- [Pixmap](#) & [GetPixmap](#) ()
- const [Pixmap](#) & [GetPixmap](#) () const
- void [SetOutputDimensions](#) (const unsigned int dims[2])  
*Set Target dimension of output Icon.*
- void [SetOutsideValuePixel](#) (double v)
- void [SetPixelMinMax](#) (double min, double max)
- void [SetPixmap](#) (const [Pixmap](#) &p)  
*Set/Get File.*

### 27.136.1 Detailed Description

[IconImageGenerator](#) This filter will generate a valid Icon from the Pixel Data element (an instance of [Pixmap](#)). To generate a valid Icon, one is only allowed the following Photometric Interpretation:

- MONOCHROME1
- MONOCHROME2
- PALETTE\_COLOR

The Pixel Bits Allocated is restricted to 8bits, therefore 16 bits image needs to be rescaled. By default the filter will use the full scalar range of 16bits image to rescale to unsigned 8bits. This may not be ideal for some situation, in which case the API `SetPixelMinMax` can be used to overwrite the default min,max interval used.

See also

[ImageReader](#)

Examples:

[ExtractIconFromFile.cxx](#).

### 27.136.2 Constructor & Destructor Documentation

27.136.2.1 `gdcm::IconImageGenerator::IconImageGenerator ( )`

27.136.2.2 `gdcm::IconImageGenerator::~~IconImageGenerator ( )`

### 27.136.3 Member Function Documentation

27.136.3.1 `void gdcm::IconImageGenerator::AutoPixelMinMax ( bool b )`

Instead of explicitly specifying the min/max value for the rescale operation, let the internal mechanism compute the min/max of icon and rescale to best appropriate.

Examples:

[ExtractIconFromFile.cxx](#).

27.136.3.2 `void gdcm::IconImageGenerator::ConvertRGBToPaletteColor ( bool b )`

Converting from RGB to PALETTE\_COLOR can be a slow operation. However DICOM standard requires that color icon be described as palette. Set this boolean to false only if you understand the consequences. default value is true, false generates invalid Icon [Image](#) Sequence

27.136.3.3 `bool gdcm::IconImageGenerator::Generate ( )`

Generate Icon.

Examples:

[ExtractIconFromFile.cxx](#).

27.136.3.4 `const IconImage& gdcm::IconImageGenerator::GetIconImage ( ) const` `[inline]`

Retrieve generated Icon.

Examples:

[ExtractIconFromFile.cxx](#).

27.136.3.5 `Pixmap& gdcm::IconImageGenerator::GetPixmap ( )` `[inline]`

27.136.3.6 `const Pixmap& gdcm::IconImageGenerator::GetPixmap ( ) const` `[inline]`

27.136.3.7 `void gdcm::IconImageGenerator::SetOutputDimensions ( const unsigned int dims[2] )`

Set Target dimension of output Icon.

Examples:

[ExtractIconFromFile.cxx](#).

27.136.3.8 `void gdcm::IconImageGenerator::SetOutsideValuePixel ( double v )`

Set a pixel value that should be discarded. This happen typically for CT image, where a pixel has been used to pad outside the image (see Pixel Padding [Value](#)). Requires `AutoPixelMinMax(true)`

27.136.3.9 `void gdcm::IconImageGenerator::SetPixelMinMax ( double min, double max )`

Override default min/max to compute best rescale for 16bits -> 8bits downscale. Typically those value can be read from the `SmallestImagePixelValue` `LargestImagePixelValue` DICOM attribute.

27.136.3.10 `void gdcm::IconImageGenerator::SetPixmap ( const Pixmap & p )` `[inline]`

Set/Get [File](#).

Examples:

[ExtractIconFromFile.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmIconImageGenerator.h](#)

## 27.137 gdcm::ignore\_char Struct Reference

```
#include <gdcmElement.h>
```

### Public Member Functions

- [ignore\\_char](#) (char c)



## Public Attributes

- char [m\\_char](#)

## 27.137.1 Constructor & Destructor Documentation

27.137.1.1 `gdcm::ignore_char::ignore_char ( char c )` `[inline]`

## 27.137.2 Member Data Documentation

27.137.2.1 `char gdcm::ignore_char::m_char`

Referenced by `gdcm::operator>>()`.

The documentation for this struct was generated from the following file:


- [gdcmElement.h](#)

## 27.138 gdcm::Image Class Reference

[Image](#) This is the container for an [Image](#) in the general sense. From this container you should be able to request information like:

```
#include <gdcmImage.h>
```

Inheritance diagram for `gdcm::Image`:

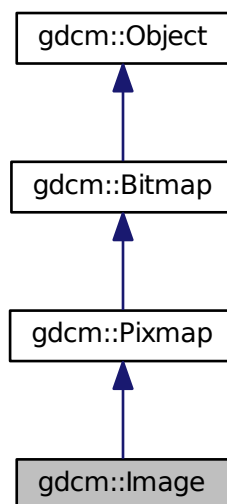



- Origin
- Dimension
- [PixelFormat](#) ... But also to retrieve the image as a raw buffer (char \*) Since we have to deal with both RAW data and JPEG stream (which internally encode all the above information) this API might seems redundant. One way to solve that would be to subclass [Image](#) with [JPEGImage](#) which would from the stream extract the header info and fill it to please [Image](#)...well except origin for instance

Basically you can see it as a storage for the Pixel Data element (7fe0,0010).

#### Warning

This class does some heuristics to guess the [Spacing](#) but is not compatible with DICOM CP-586. In case of doubt use [PixmapReader](#) instead

#### See also

[ImageReader](#) [PixmapReader](#)

#### Examples:

[CompressImage.cxx](#), [ConvertToQImage.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [ExtractIconFromFile.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenFakeImage.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), [iU22tomultisc.cxx](#), [PatchFile.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

## 27.138.2 Constructor & Destructor Documentation

27.138.2.1 `gdcm::Image::Image ( )` `[inline]`

27.138.2.2 `gdcm::Image::~~Image ( )` `[inline]`

## 27.138.3 Member Function Documentation

27.138.3.1 `const double* gdcm::Image::GetDirectionCosines ( )` `const`

Return a 6-tuples specifying the direction cosines A default value of (1,0,0,0,1,0) will be return when the direction cosines was not specified.

27.138.3.2 `double gdcm::Image::GetDirectionCosines ( unsigned int idx )` `const`

27.138.3.3 `double gdcm::Image::GetIntercept ( )` `const` `[inline]`

27.138.3.4 `const double* gdcm::Image::GetOrigin ( )` `const`

Return a 3-tuples specifying the origin Will return (0,0,0) if the origin was not specified.

#### Examples:

[HelloVizWorld.cxx](#).

27.138.3.5 `double gdc::Image::GetOrigin ( unsigned int idx ) const`

27.138.3.6 `double gdc::Image::GetSlope ( ) const` `[inline]`

27.138.3.7 `const double* gdc::Image::GetSpacing ( ) const`

Return a 3-tuples specifying the spacing NOTE: 3rd value can be an arbitrary 1 value when the spacing was not specified (ex. 2D image). WARNING: when the spacing is not specifier, a default value of 1 will be returned

27.138.3.8 `double gdc::Image::GetSpacing ( unsigned int idx ) const`

27.138.3.9 `void gdc::Image::Print ( std::ostream & os ) const` `[virtual]`

print

Reimplemented from [gdc::Bitmap](#).

Examples:

[CompressImage.cxx](#), and [PatchFile.cxx](#).

27.138.3.10 `void gdc::Image::SetDirectionCosines ( const float * dircos )`

27.138.3.11 `void gdc::Image::SetDirectionCosines ( const double * dircos )`

27.138.3.12 `void gdc::Image::SetDirectionCosines ( unsigned int idx, double dircos )`

27.138.3.13 `void gdc::Image::SetIntercept ( double intercept )` `[inline]`

intercept

27.138.3.14 `void gdc::Image::SetOrigin ( const float * ori )`

27.138.3.15 `void gdc::Image::SetOrigin ( const double * ori )`

27.138.3.16 `void gdc::Image::SetOrigin ( unsigned int idx, double ori )`

27.138.3.17 `void gdc::Image::SetSlope ( double slope )` `[inline]`

slope

27.138.3.18 `void gdc::Image::SetSpacing ( const double * spacing )`

Examples:

[csa2img.cxx](#), and [iU22tomultisc.cxx](#).

27.138.3.19 void gdcm::Image::SetSpacing ( unsigned int *idx*, double *spacing* )

The documentation for this class was generated from the following file:


- [gdcmImage.h](#)

## 27.139 gdcm::ImageApplyLookupTable Class Reference


[ImageApplyLookupTable](#) class It applies the LUT the PixelData (only PALETTE\_COLOR images) Output will be a [PhotometricInterpretation](#)=RGB image.

```
#include <gdcmImageApplyLookupTable.h>
```

Inheritance diagram for gdcm::ImageApplyLookupTable:

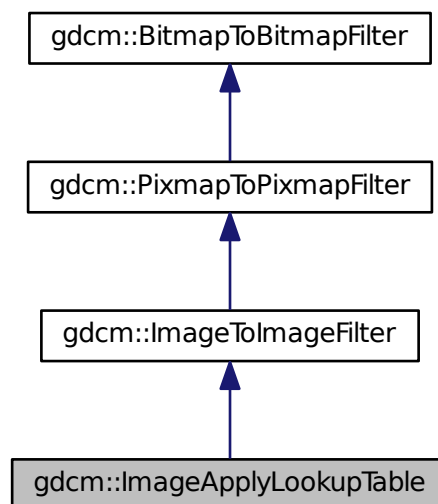

Collaboration diagram for `gdcm::ImageApplyLookupTable`:

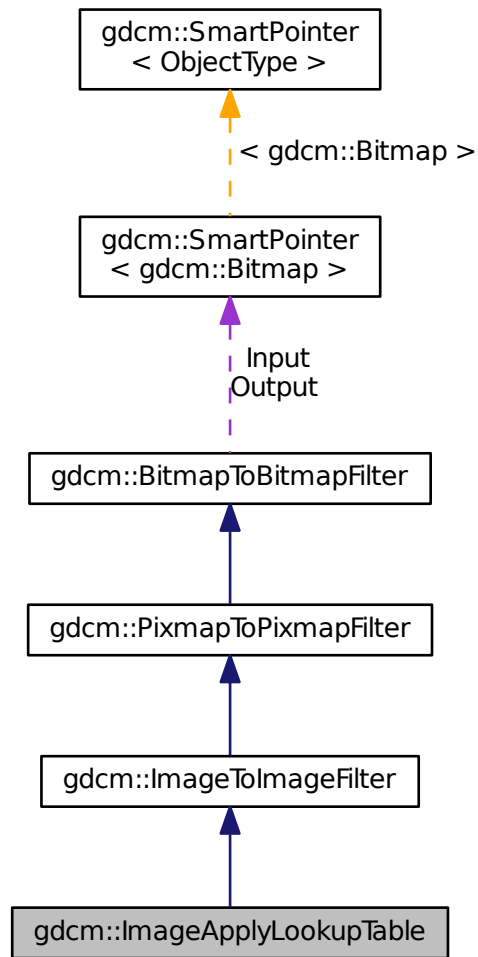

### Public Member Functions

- [ImageApplyLookupTable \(\)](#)
- [~ImageApplyLookupTable \(\)](#)
- `bool` [Apply \(\)](#)

*Apply.*

### Additional Inherited Members

### 27.139.1 Detailed Description

[ImageApplyLookupTable](#) class It applies the LUT the PixelData (only PALETTE\_COLOR images) Output will be a [PhotometricInterpretation](#)=RGB image.

### 27.139.2 Constructor & Destructor Documentation

27.139.2.1 `gdcm::ImageApplyLookupTable::ImageApplyLookupTable ( )` `[inline]`

27.139.2.2 `gdcm::ImageApplyLookupTable::~~ImageApplyLookupTable ( )` `[inline]`

### 27.139.3 Member Function Documentation

27.139.3.1 `bool gdcm::ImageApplyLookupTable::Apply ( )`

Apply.

The documentation for this class was generated from the following file:


- [gdcmImageApplyLookupTable.h](#)

## 27.140 gdcm::ImageChangePhotometricInterpretation Class Reference


[ImageChangePhotometricInterpretation](#) class Class to change the Photometric Interpretation of an input DICOM.

```
#include <gdcmImageChangePhotometricInterpretation.h>
```

Inheritance diagram for `gdcm::ImageChangePhotometricInterpretation`:

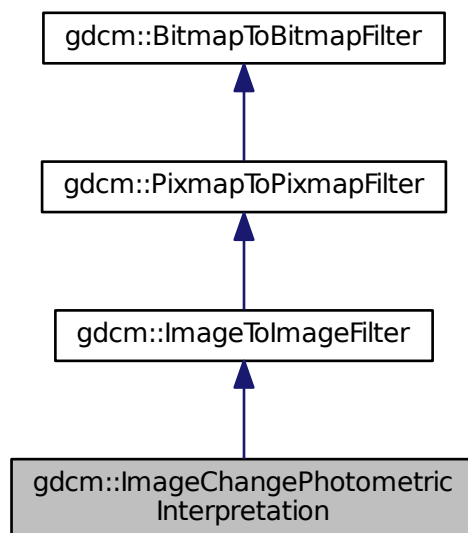

Collaboration diagram for `gdcm::ImageChangePhotometricInterpretation`:

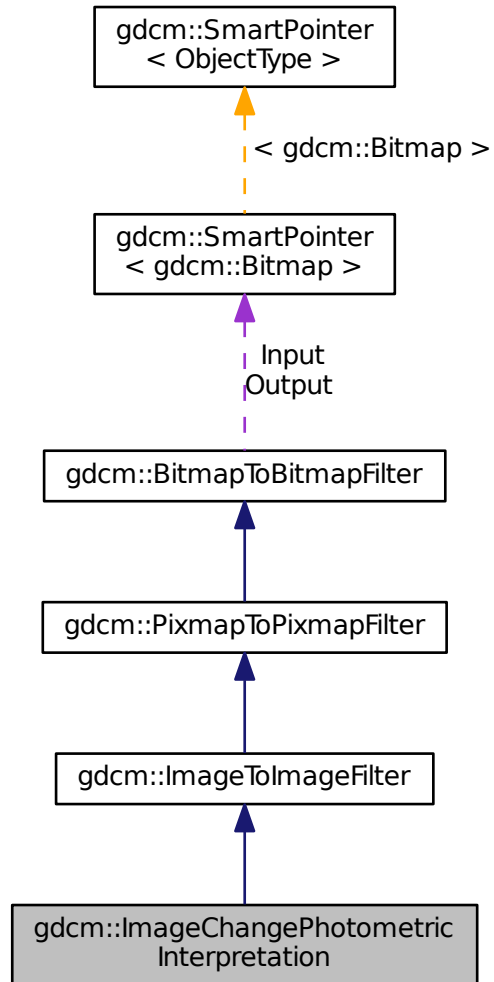

## Public Member Functions

- [ImageChangePhotometricInterpretation \(\)](#)
- [~ImageChangePhotometricInterpretation \(\)](#)
- [bool Change \(\)](#)  
*Change.*
- [const PhotometricInterpretation & GetPhotometricInterpretation \(\) const](#)
- [void SetPhotometricInterpretation \(PhotometricInterpretation const &pi\)](#)  
*Set/Get requested PhotometricInterpretation.*



## Static Public Member Functions

- template<typename T >  
static void [RGB2YBR](#) (T ybr[3], const T rgb[3])  
*colorspace conversion (based on CCIR Recommendation 601-2)*
- template<typename T >  
static void [YBR2RGB](#) (T rgb[3], const T ybr[3])

## Protected Member Functions

- bool [ChangeMonochrome](#) ()

## Additional Inherited Members

### 27.140.1 Detailed Description

[ImageChangePhotometricInterpretation](#) class Class to change the Photometric Interpretation of an input DICOM.

### 27.140.2 Constructor & Destructor Documentation

27.140.2.1 `gdcm::ImageChangePhotometricInterpretation::ImageChangePhotometricInterpretation ( )` `[inline]`

27.140.2.2 `gdcm::ImageChangePhotometricInterpretation::~~ImageChangePhotometricInterpretation ( )` `[inline]`

### 27.140.3 Member Function Documentation

27.140.3.1 `bool gdcm::ImageChangePhotometricInterpretation::Change ( )`

Change.

27.140.3.2 `bool gdcm::ImageChangePhotometricInterpretation::ChangeMonochrome ( )` `[protected]`

27.140.3.3 `const PhotometricInterpretation& gdcm::ImageChangePhotometricInterpretation::GetPhotometricInterpretation ( )`  
`const` `[inline]`

27.140.3.4 `template<typename T > void gdcm::ImageChangePhotometricInterpretation::RGB2YBR ( T ybr[3], const T rgb[3] )`  
`[static]`

colorspace conversion (based on CCIR Recommendation 601-2)

27.140.3.5 `void gdcm::ImageChangePhotometricInterpretation::SetPhotometricInterpretation ( PhotometricInterpretation const & pi )` `[inline]`

Set/Get requested [PhotometricInterpretation](#).

27.140.3.6 `template<typename T > void gdcmm::ImageChangePhotometricInterpretation::YBR2RGB ( T rgb[3], const T ybr[3] )`  
`[static]`

The documentation for this class was generated from the following file:


- [gdcmImageChangePhotometricInterpretation.h](#)

## 27.141 gdcm::ImageChangePlanarConfiguration Class Reference


[ImageChangePlanarConfiguration](#) class Class to change the Planar configuration of an input DICOM By default it will change into the more usual representation: PlanarConfiguration = 0.

```
#include <gdcmImageChangePlanarConfiguration.h>
```

Inheritance diagram for `gdcm::ImageChangePlanarConfiguration`:

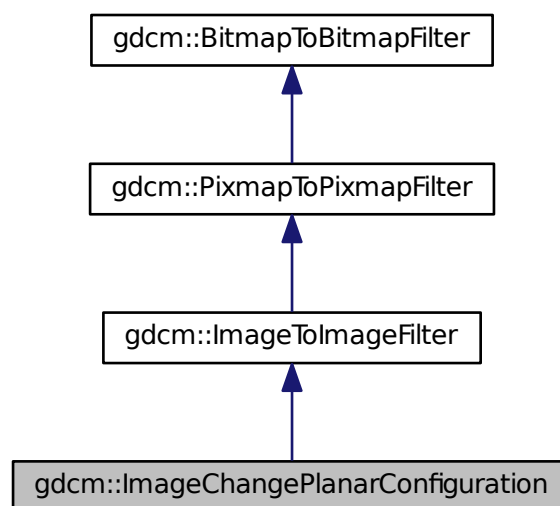

Collaboration diagram for gdcm::ImageChangePlanarConfiguration:

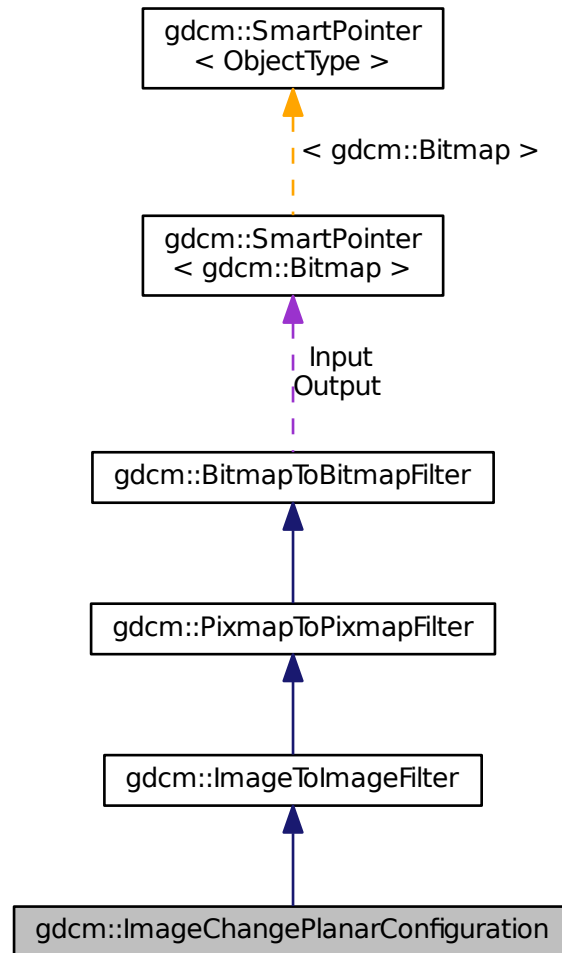

## Public Member Functions

- [ImageChangePlanarConfiguration](#) ()
- [~ImageChangePlanarConfiguration](#) ()
- [bool Change](#) ()  
*Change.*
- [unsigned int GetPlanarConfiguration](#) () const
- [void SetPlanarConfiguration](#) (unsigned int pc)  
*Set/Get requested PlanarConfiguration.*

## Static Public Member Functions

- `template<typename T >`  
`static size_t RGBPixelsToRGBPlanes (T *r, T *g, T *b, const T *rgb, size_t s)`
- `template<typename T >`  
`static size_t RGBPlanesToRGBPixels (T *out, const T *r, const T *g, const T *b, size_t s)`

## Additional Inherited Members

### 27.141.1 Detailed Description

[ImageChangePlanarConfiguration](#) class Class to change the Planar configuration of an input DICOM By default it will change into the more usual representation: PlanarConfiguration = 0.

### 27.141.2 Constructor & Destructor Documentation

27.141.2.1 `gdcm::ImageChangePlanarConfiguration::ImageChangePlanarConfiguration ( )` `[inline]`

27.141.2.2 `gdcm::ImageChangePlanarConfiguration::~~ImageChangePlanarConfiguration ( )` `[inline]`

### 27.141.3 Member Function Documentation

27.141.3.1 `bool gdcm::ImageChangePlanarConfiguration::Change ( )`

Change.

27.141.3.2 `unsigned int gdcm::ImageChangePlanarConfiguration::GetPlanarConfiguration ( ) const` `[inline]`

27.141.3.3 `template<typename T > size_t gdcm::ImageChangePlanarConfiguration::RGBPixelsToRGBPlanes ( T * r, T * g, T * b, const T * rgb, size_t s )` `[static]`

Convert a regular RGB pixel image (R,G,B,R,G,B...) into a planar R,G,B image (R,R...,G,G...B,B)

#### Warning

this works on a frame basis, you need to loop over all frames in multiple frames image to apply this function

27.141.3.4 `template<typename T > size_t gdcm::ImageChangePlanarConfiguration::RGBPlanesToRGBPixels ( T * out, const T * r, const T * g, const T * b, size_t s )` `[static]`

s is the size of one plane (r,g or b). Thus the output buffer needs to be at least 3\*s bytes long s can be seen as the number of RGB pixels in the output

27.141.3.5 `void gdcm::ImageChangePlanarConfiguration::SetPlanarConfiguration ( unsigned int pc )` `[inline]`

Set/Get requested PlanarConfiguration.

The documentation for this class was generated from the following file:


- [gdcmImageChangePlanarConfiguration.h](#)

## 27.142 gdcm::ImageChangeTransferSyntax Class Reference


[ImageChangeTransferSyntax](#) class Class to change the transfer syntax of an input DICOM.

```
#include <gdcmImageChangeTransferSyntax.h>
```

Inheritance diagram for gdcm::ImageChangeTransferSyntax:

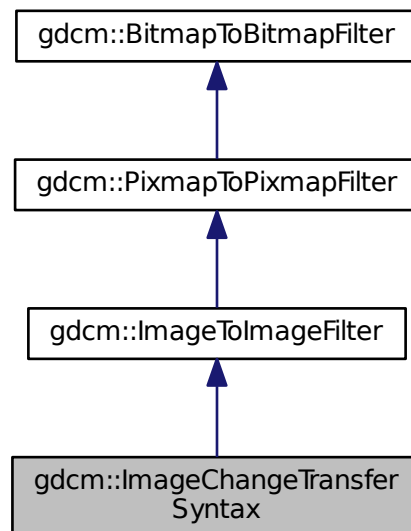

Collaboration diagram for `gdcm::ImageChangeTransferSyntax`:

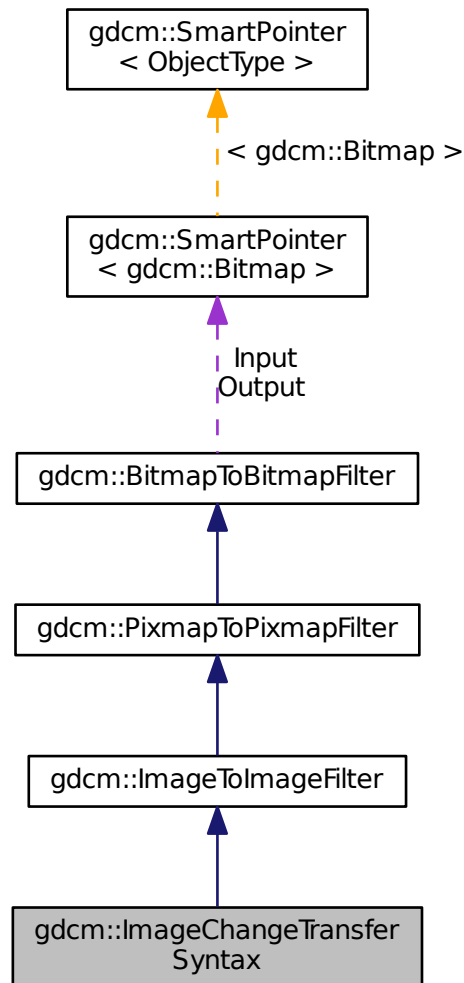

## Public Member Functions

- [ImageChangeTransferSyntax](#) ()
- [~ImageChangeTransferSyntax](#) ()
- [bool Change](#) ()
  - Change.*
- [const TransferSyntax & GetTransferSyntax](#) () const
  - Get Transfer Syntax.*
- [void SetCompressIconImage](#) (bool b)
- [void SetForce](#) (bool f)
- [void SetTransferSyntax](#) (const [TransferSyntax](#) &ts)

*Set target Transfer Syntax.*

- void [SetUserCodec](#) ([ImageCodec](#) \*ic)

## Protected Member Functions

- bool [TryJPEG2000Codec](#) (const [DataElement](#) &pixelde, [Bitmap](#) const &input, [Bitmap](#) &output)
- bool [TryJPEGCodec](#) (const [DataElement](#) &pixelde, [Bitmap](#) const &input, [Bitmap](#) &output)
- bool [TryJPEGLSCoec](#) (const [DataElement](#) &pixelde, [Bitmap](#) const &input, [Bitmap](#) &output)
- bool [TryRAWCodec](#) (const [DataElement](#) &pixelde, [Bitmap](#) const &input, [Bitmap](#) &output)
- bool [TryRLECodec](#) (const [DataElement](#) &pixelde, [Bitmap](#) const &input, [Bitmap](#) &output)

## Additional Inherited Members

### 27.142.1 Detailed Description

[ImageChangeTransferSyntax](#) class Class to change the transfer syntax of an input DICOM.

If only Force param is set but no input [TransferSyntax](#) is set, it is assumed that user only wants to inspect encapsulated stream (advanced dev. option).

When using UserCodec it is very important that the [TransferSyntax](#) (as set in [SetTransferSyntax](#)) is actually understood by UserCodec (ie. `UserCodec->CanCode( TransferSyntax )` ). Otherwise the behavior is to use a default codec.

See also

[JPEGCodec](#) [JPEGLSCoec](#) [JPEG2000Codec](#)

Examples:

[CompressImage.cxx](#).

### 27.142.2 Constructor & Destructor Documentation

27.142.2.1 `gdcm::ImageChangeTransferSyntax::ImageChangeTransferSyntax ( )` `[inline]`

27.142.2.2 `gdcm::ImageChangeTransferSyntax::~ImageChangeTransferSyntax ( )` `[inline]`

### 27.142.3 Member Function Documentation

27.142.3.1 `bool gdcm::ImageChangeTransferSyntax::Change ( )`

Change.

Examples:

[CompressImage.cxx](#).

27.142.3.2 `const TransferSyntax& gdcm::ImageChangeTransferSyntax::GetTransferSyntax ( ) const` `[inline]`

Get Transfer Syntax.

27.142.3.3 `void gdcmm::ImageChangeTransferSyntax::SetCompressIconImage ( bool b )` `[inline]`

Decide whether or not to also compress the Icon [Image](#) using the same Transfer Syntax. Default is to simply decompress icon image

27.142.3.4 `void gdcmm::ImageChangeTransferSyntax::SetForce ( bool f )` `[inline]`

When target Transfer Syntax is identical to input target syntax, no operation is actually done. This is an issue when someone wants to re-compress using GDCM internal implementation a JPEG (for example) image

27.142.3.5 `void gdcmm::ImageChangeTransferSyntax::SetTransferSyntax ( const TransferSyntax & ts )` `[inline]`

Set target Transfer Syntax.

Examples:

[CompressImage.cxx](#).

27.142.3.6 `void gdcmm::ImageChangeTransferSyntax::SetUserCodec ( ImageCodec * ic )` `[inline]`

Allow user to specify exactly which codec to use. this is needed to specify special qualities or compression option.

Warning

if the codec '*ic*' is not compatible with the [TransferSyntax](#) requested, it will not be used. It is the user responsibility to check that `UserCodec->CanCode( TransferSyntax )`

27.142.3.7 `bool gdcmm::ImageChangeTransferSyntax::TryJPEG2000Codec ( const DataElement & pixelde, Bitmap const & input, Bitmap & output )` `[protected]`

27.142.3.8 `bool gdcmm::ImageChangeTransferSyntax::TryJPEGCodec ( const DataElement & pixelde, Bitmap const & input, Bitmap & output )` `[protected]`

27.142.3.9 `bool gdcmm::ImageChangeTransferSyntax::TryJPEGLSCodec ( const DataElement & pixelde, Bitmap const & input, Bitmap & output )` `[protected]`

27.142.3.10 `bool gdcmm::ImageChangeTransferSyntax::TryRAWCodec ( const DataElement & pixelde, Bitmap const & input, Bitmap & output )` `[protected]`

27.142.3.11 `bool gdcmm::ImageChangeTransferSyntax::TryRLECodec ( const DataElement & pixelde, Bitmap const & input, Bitmap & output )` `[protected]`

The documentation for this class was generated from the following file:

- [gdcmmImageChangeTransferSyntax.h](#)




## 27.143 gdcm::ImageCodec Class Reference


[ImageCodec](#).

```
#include <gdcmImageCodec.h>
```

Inheritance diagram for gdcm::ImageCodec:

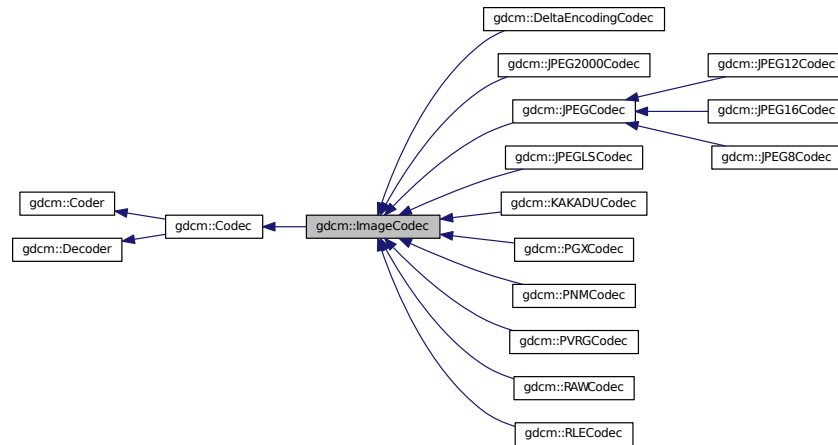

Collaboration diagram for gdcm::ImageCodec:

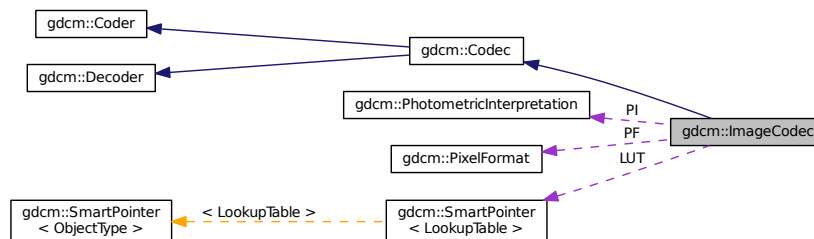

### Public Member Functions

- [ImageCodec](#) ()
- [~ImageCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &) const  
Return whether this coder support this transfer syntax (can code it)
- bool [CanDecode](#) ([TransferSyntax](#) const &) const  
Return whether this decoder support this transfer syntax (can decode it)
- virtual [ImageCodec](#) \* [Clone](#) () const =0
- bool [Decode](#) ([DataElement](#) const &is\_, [DataElement](#) &os)  
Decode.

- const unsigned int \* [GetDimensions](#) () const
- virtual bool [GetHeaderInfo](#) (std::istream &is\_, [TransferSyntax](#) &ts)
- bool [GetLossyFlag](#) () const
- const [LookupTable](#) & [GetLUT](#) () const
- bool [GetNeedByteSwap](#) () const
- unsigned int [GetNumberOfDimensions](#) () const
- const [PhotometricInterpretation](#) & [GetPhotometricInterpretation](#) () const
- [PixelFormat](#) & [GetPixelFormat](#) ()
- const [PixelFormat](#) & [GetPixelFormat](#) () const
- unsigned int [GetPlanarConfiguration](#) () const
- bool [IsLossy](#) () const
- void [SetDimensions](#) (const unsigned int d[3])
- void [SetDimensions](#) (const std::vector< unsigned int > &d)
- void [SetLossyFlag](#) (bool l)
- void [SetLUT](#) ([LookupTable](#) const &lut)
- void [SetNeedByteSwap](#) (bool b)
- void [SetNeedOverlayCleanup](#) (bool b)
- void [SetNumberOfDimensions](#) (unsigned int dim)
- void [SetPhotometricInterpretation](#) ([PhotometricInterpretation](#) const &pi)
- virtual void [SetPixelFormat](#) ([PixelFormat](#) const &pf)
- void [SetPlanarConfiguration](#) (unsigned int pc)

## Protected Types

- typedef [SmartPointer](#)< [LookupTable](#) > [LUTPtr](#)

## Protected Member Functions

- virtual bool [AppendFrameEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- virtual bool [AppendRowEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [DecodeByStreams](#) (std::istream &is\_, std::ostream &os)
- bool [DoByteSwap](#) (std::istream &is\_, std::ostream &os)
- bool [DoInvertMonochrome](#) (std::istream &is\_, std::ostream &os)
- bool [DoOverlayCleanup](#) (std::istream &is\_, std::ostream &os)
- bool [DoPaddedCompositePixelCode](#) (std::istream &is\_, std::ostream &os)
- bool [DoPlanarConfiguration](#) (std::istream &is\_, std::ostream &os)
- bool [DoSimpleCopy](#) (std::istream &is\_, std::ostream &os)
- bool [DoYBR](#) (std::istream &is\_, std::ostream &os)
- virtual bool [IsFrameEncoder](#) ()
- virtual bool [IsRowEncoder](#) ()
- virtual bool [IsValid](#) ([PhotometricInterpretation](#) const &pi)
- virtual bool [StartEncode](#) (std::ostream &os)
- virtual bool [StopEncode](#) (std::ostream &os)

## Protected Attributes

- unsigned int [Dimensions](#) [3]
- bool [LossyFlag](#)
- [LUTPtr](#) LUT
- bool [NeedByteSwap](#)
- bool [NeedOverlayCleanup](#)
- unsigned int [NumberOfDimensions](#)
- [PixelFormat](#) PF
- [PhotometricInterpretation](#) PI
- unsigned int [PlanarConfiguration](#)
- bool [RequestPaddedCompositePixelCode](#)
- bool [RequestPlanarConfiguration](#)

## Friends

- class [FileChangeTransferSyntax](#)
- class [ImageChangePhotometricInterpretation](#)

### 27.143.1 Detailed Description

[ImageCodec](#).

#### Note

Main codec, this is a central place for all implementation

### 27.143.2 Member Typedef Documentation

27.143.2.1 `typedef SmartPointer<LookupTable> gdcm::ImageCodec::LUTPtr` [protected]

### 27.143.3 Constructor & Destructor Documentation

27.143.3.1 `gdcm::ImageCodec::ImageCodec ( )`

27.143.3.2 `gdcm::ImageCodec::~~ImageCodec ( )`

### 27.143.4 Member Function Documentation

27.143.4.1 `virtual bool gdcm::ImageCodec::AppendFrameEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

27.143.4.2 `virtual bool gdcm::ImageCodec::AppendRowEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

27.143.4.3 `bool gdcm::ImageCodec::CanCode ( TransferSyntax const & ) const` `[inline],[virtual]`

Return whether this coder support this transfer syntax (can code it)

Implements [gdcm::Coder](#).

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::JPEGGLSCodec](#), [gdcm::PNMCodec](#), [gdcm::PGXCodec](#), [gdcm::KAKADUCodec](#), and [gdcm::RAWCodec](#).

27.143.4.4 `bool gdcm::ImageCodec::CanDecode ( TransferSyntax const & ) const` `[inline],[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Implements [gdcm::Decoder](#).

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::JPEGGLSCodec](#), [gdcm::PNMCodec](#), [gdcm::RAWCodec](#), [gdcm::PGXCodec](#), and [gdcm::KAKADUCodec](#).

27.143.4.5 `virtual ImageCodec* gdcm::ImageCodec::Clone ( ) const` `[pure virtual]`

Implemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::JPEGGLSCodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::PNMCodec](#), [gdcm::RAWCodec](#), [gdcm::KAKADUCodec](#), and [gdcm::PGXCodec](#).

27.143.4.6 `bool gdcm::ImageCodec::Decode ( DataElement const & , DataElement & )` `[virtual]`

Decode.

Reimplemented from [gdcm::Decoder](#).

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::JPEGGLSCodec](#), [gdcm::PVRGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::KAKADUCodec](#), and [gdcm::RAWCodec](#).

27.143.4.7 `bool gdcm::ImageCodec::DecodeByStreams ( std::istream & is_ , std::ostream & os )` `[protected],[virtual]`

Reimplemented from [gdcm::Decoder](#).

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::RLECodec](#), [gdcm::RAWCodec](#), [gdcm::JPEGGLSCodec](#), [gdcm::JPEG16Codec](#), and [gdcm::JPEG8Codec](#).

27.143.4.8 `bool gdcm::ImageCodec::DoByteSwap ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.9 `bool gdcm::ImageCodec::DoInvertMonochrome ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.10 `bool gdcm::ImageCodec::DoOverlayCleanup ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.11 `bool gdcm::ImageCodec::DoPaddedCompositePixelCode ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.12 `bool gdcm::ImageCodec::DoPlanarConfiguration ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.13 `bool gdcm::ImageCodec::DoSimpleCopy ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.14 `bool gdcm::ImageCodec::DoYBR ( std::istream & is_ , std::ostream & os )` `[protected]`

27.143.4.15 `const unsigned int* gdcm::ImageCodec::GetDimensions ( ) const` [inline]

27.143.4.16 `virtual bool gdcm::ImageCodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::RLECodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), [gdcm::PNMCodec](#), [gdcm::JPEG12Codec](#), [gdcm::JPEG16Codec](#), [gdcm::JPEG8Codec](#), [gdcm::RAWCodec](#), and [gdcm::PGXCodec](#).

27.143.4.17 `bool gdcm::ImageCodec::GetLossyFlag ( ) const`

27.143.4.18 `const LookupTable& gdcm::ImageCodec::GetLUT ( ) const` [inline]

27.143.4.19 `bool gdcm::ImageCodec::GetNeedByteSwap ( ) const` [inline]

27.143.4.20 `unsigned int gdcm::ImageCodec::GetNumberOfDimensions ( ) const`

27.143.4.21 `const PhotometricInterpretation& gdcm::ImageCodec::GetPhotometricInterpretation ( ) const`

27.143.4.22 `PixelFormat& gdcm::ImageCodec::GetPixelFormat ( )` [inline]

Examples:

[GetJPEGSamplePrecision.cxx](#).

27.143.4.23 `const PixelFormat& gdcm::ImageCodec::GetPixelFormat ( ) const` [inline]

27.143.4.24 `unsigned int gdcm::ImageCodec::GetPlanarConfiguration ( ) const` [inline]

27.143.4.25 `virtual bool gdcm::ImageCodec::IsFrameEncoder ( )` [protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

27.143.4.26 `bool gdcm::ImageCodec::IsLossy ( ) const`

27.143.4.27 `virtual bool gdcm::ImageCodec::IsRowEncoder ( )` [protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

27.143.4.28 `virtual bool gdcm::ImageCodec::IsValid ( PhotometricInterpretation const & pi )` [protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#).

27.143.4.29 `void gdcm::ImageCodec::SetDimensions ( const unsigned int d[3] )`

Examples:

[ExtractIconFromFile.cxx](#).

27.143.4.30 void `gdcm::ImageCodec::SetDimensions` ( const `std::vector< unsigned int >` & *d* )

27.143.4.31 void `gdcm::ImageCodec::SetLossyFlag` ( bool *l* )

27.143.4.32 void `gdcm::ImageCodec::SetLUT` ( `LookupTable` const & *lut* ) [inline]

Examples:

[ExtractIconFromFile.cxx](#).

27.143.4.33 void `gdcm::ImageCodec::SetNeedByteSwap` ( bool *b* ) [inline]

27.143.4.34 void `gdcm::ImageCodec::SetNeedOverlayCleanup` ( bool *b* ) [inline]

27.143.4.35 void `gdcm::ImageCodec::SetNumberOfDimensions` ( unsigned int *dim* )

27.143.4.36 void `gdcm::ImageCodec::SetPhotometricInterpretation` ( `PhotometricInterpretation` const & *pi* )

Examples:

[ExtractIconFromFile.cxx](#).

27.143.4.37 virtual void `gdcm::ImageCodec::SetPixelFormat` ( `PixelFormat` const & *pf* ) [inline],[virtual]

Reimplemented in [gdcm::JPEGCodec](#).

Examples:

[ExtractIconFromFile.cxx](#).

27.143.4.38 void `gdcm::ImageCodec::SetPlanarConfiguration` ( unsigned int *pc* ) [inline]

27.143.4.39 virtual bool `gdcm::ImageCodec::StartEncode` ( `std::ostream` & *os* ) [protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

27.143.4.40 virtual bool `gdcm::ImageCodec::StopEncode` ( `std::ostream` & *os* ) [protected],[virtual]

Reimplemented in [gdcm::JPEGCodec](#), [gdcm::JPEGLSCodec](#), [gdcm::JPEG2000Codec](#), and [gdcm::RLECodec](#).

## 27.143.5 Friends And Related Function Documentation

27.143.5.1 friend class `FileChangeTransferSyntax` [friend]

This is a high level API to encode in a streaming fashion. Each plugin will handle differently the caching mechanism so that a limited memory is used when compressing dataset. [Codec](#) will fall into two categories:

- Full row encoder: only a single scanline (row) of data is needed to be loaded at a time;
- Full frame encoder (default): a complete frame (row x col) is needed to be loaded at a time

27.143.5.2 friend class ImageChangePhotometricInterpretation [friend]

### 27.143.6 Member Data Documentation

27.143.6.1 unsigned int gdcm::ImageCodec::Dimensions[3] [protected]

27.143.6.2 bool gdcm::ImageCodec::LossyFlag [protected]

27.143.6.3 LUTPtr gdcm::ImageCodec::LUT [protected]

27.143.6.4 bool gdcm::ImageCodec::NeedByteSwap [protected]

27.143.6.5 bool gdcm::ImageCodec::NeedOverlayCleanup [protected]

27.143.6.6 unsigned int gdcm::ImageCodec::NumberOfDimensions [protected]

27.143.6.7 PixelFormat gdcm::ImageCodec::PF [protected]

27.143.6.8 PhotometricInterpretation gdcm::ImageCodec::PI [protected]

27.143.6.9 unsigned int gdcm::ImageCodec::PlanarConfiguration [protected]

27.143.6.10 bool gdcm::ImageCodec::RequestPaddedCompositePixelCode [protected]

27.143.6.11 bool gdcm::ImageCodec::RequestPlanarConfiguration [protected]

The documentation for this class was generated from the following file:

- [gdcmImageCodec.h](#)

## 27.144 gdcm::ImageConverter Class Reference

[Image](#) Converter.

```
#include <gdcmImageConverter.h>
```

### Public Member Functions

- [ImageConverter](#) ()
- [~ImageConverter](#) ()
- void [Convert](#) ()
- const [Image](#) & [GetOutput](#) () const
- void [SetInput](#) ([Image](#) const &input)

### 27.144.1 Detailed Description

[Image](#) Converter.

**Note**

This is the class used to convert from one [Image](#) to another. This is typically used to convert let say YBR JPEG compressed [Image](#) to a RAW RGB [Image](#). So that the buffer can be directly pass to third party application. This filter is application level and not integrated directly in GDCM

**27.144.2 Constructor & Destructor Documentation**

27.144.2.1 `gdcm::ImageConverter::ImageConverter ( )`

27.144.2.2 `gdcm::ImageConverter::~~ImageConverter ( )`

**27.144.3 Member Function Documentation**

27.144.3.1 `void gdcm::ImageConverter::Convert ( )`

27.144.3.2 `const Image& gdcm::ImageConverter::GetOutput ( ) const`

27.144.3.3 `void gdcm::ImageConverter::SetInput ( Image const & input )`

The documentation for this class was generated from the following file:

- [gdcmImageConverter.h](#)


**27.145 gdcm::ImageFragmentSplitter Class Reference**

[ImageFragmentSplitter](#) class For single frame image, DICOM standard allow splitting the frame into multiple fragments.


```
#include <gdcmImageFragmentSplitter.h>
```



Inheritance diagram for gdcm::ImageFragmentSplitter:

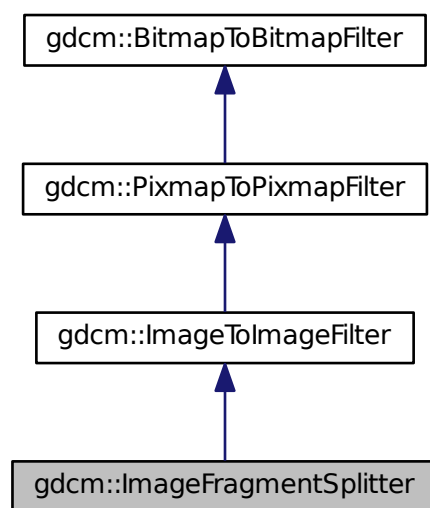

Collaboration diagram for `gdcm::ImageFragmentSplitter`:

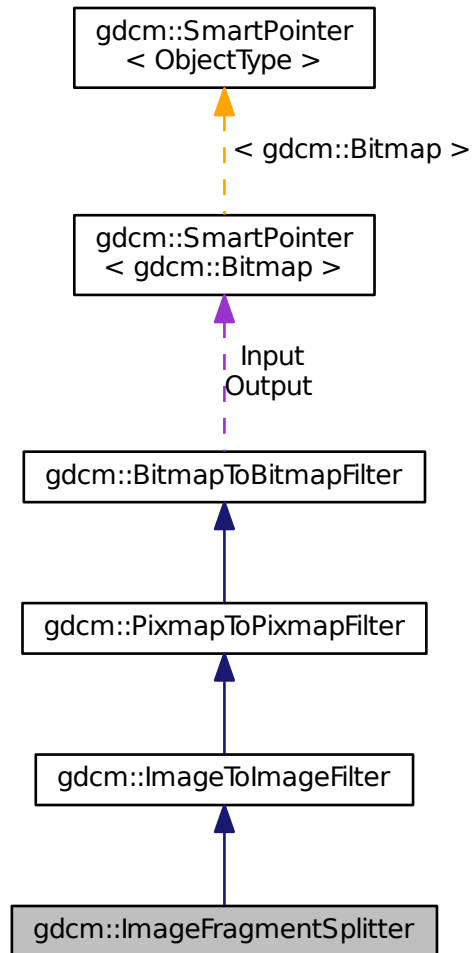

## Public Member Functions

- [ImageFragmentSplitter](#) ()
- [~ImageFragmentSplitter](#) ()
- unsigned int [GetFragmentSizeMax](#) () const
- void [SetForce](#) (bool f)
- void [SetFragmentSizeMax](#) (unsigned int fragsize)

*FragmentSizeMax needs to be an even number.*

- bool [Split](#) ()

*Split.*

## Additional Inherited Members

### 27.145.1 Detailed Description

[ImageFragmentSplitter](#) class For single frame image, DICOM standard allow splitting the frame into multiple fragments.

### 27.145.2 Constructor & Destructor Documentation

27.145.2.1 `gdcm::ImageFragmentSplitter::ImageFragmentSplitter ( )` [\[inline\]](#)

27.145.2.2 `gdcm::ImageFragmentSplitter::~~ImageFragmentSplitter ( )` [\[inline\]](#)

### 27.145.3 Member Function Documentation

27.145.3.1 `unsigned int gdcm::ImageFragmentSplitter::GetFragmentSizeMax ( ) const` [\[inline\]](#)

27.145.3.2 `void gdcm::ImageFragmentSplitter::SetForce ( bool f )` [\[inline\]](#)

When file already has all it's segment < FragmentSizeMax there is not need to run the filter. Unless the user explicitly say 'force' recomputation !

27.145.3.3 `void gdcm::ImageFragmentSplitter::SetFragmentSizeMax ( unsigned int fragsize )`

FragmentSizeMax needs to be an even number.

27.145.3.4 `bool gdcm::ImageFragmentSplitter::Split ( )`

Split.

The documentation for this class was generated from the following file:

- [gdcmImageFragmentSplitter.h](#)

## 27.146 gdcm::ImageHelper Class Reference

[ImageHelper](#) (internal class, not intended for user level)

```
#include <gdcmImageHelper.h>
```

### Static Public Member Functions

- static [MediaStorage](#) [ComputeMediaStorageFromModality](#) (const char \*modality, unsigned int dimension=2, [PixelFormat](#) const &pf=[PixelFormat](#)(), [PhotometricInterpretation](#) const &pi=[PhotometricInterpretation](#)(), double rescaleintercept=0, double rescaleslope=1)  
*Moved from [MediaStorage](#) here, since we need extra info stored in [PixelFormat](#) & [PhotometricInterpretation](#).*
- static bool [ComputeSpacingFromImagePositionPatient](#) (const std::vector< double > &imageposition, std::vector< double > &spacing)  
*DO NOT USE.*
- static std::vector< unsigned int > [GetDimensionsValue](#) (const [File](#) &f)

- static bool [GetDirectionCosinesFromDataSet](#) ([DataSet](#) const &ds, std::vector< double > &dircos)
- static std::vector< double > [GetDirectionCosinesValue](#) ([File](#) const &f)
- static bool [GetForcePixelSpacing](#) ()
- static bool [GetForceRescaleInterceptSlope](#) ()
- static [SmartPointer](#)< [LookupTable](#) > [GetLUT](#) ([File](#) const &f)
- static std::vector< double > [GetOriginValue](#) ([File](#) const &f)

*Set/Get Origin (IPP) from/to a file.*

- static [PhotometricInterpretation](#) [GetPhotometricInterpretationValue](#) ([File](#) const &f)
- static [PixelFormat](#) [GetPixelFormatValue](#) (const [File](#) &f)
- static unsigned int [GetPlanarConfigurationValue](#) (const [File](#) &f)
- static const [ByteValue](#) \* [GetPointerFromElement](#) ([Tag](#) const &tag, [File](#) const &f)

*Moved from [PixampReader](#) to here. Generally used for photometric interpretation.*

- static std::vector< double > [GetRescaleInterceptSlopeValue](#) ([File](#) const &f)
- static std::vector< double > [GetSpacingValue](#) ([File](#) const &f)

*Set/Get [Spacing](#) from/to a [File](#).*

- static void [SetDimensionsValue](#) ([File](#) &f, const [Image](#) &img)
- static void [SetDirectionCosinesValue](#) ([DataSet](#) &ds, const std::vector< double > &dircos)
- static void [SetForcePixelSpacing](#) (bool)
- static void [SetForceRescaleInterceptSlope](#) (bool)
- static void [SetOriginValue](#) ([DataSet](#) &ds, const [Image](#) &img)
- static void [SetRescaleInterceptSlopeValue](#) ([File](#) &f, const [Image](#) &img)
- static void [SetSpacingValue](#) ([DataSet](#) &ds, const std::vector< double > &spacing)

## Static Protected Member Functions

- static [Tag](#) [GetSpacingTagFromMediaStorage](#) ([MediaStorage](#) const &ms)
- static [Tag](#) [GetZSpacingTagFromMediaStorage](#) ([MediaStorage](#) const &ms)

### 27.146.1 Detailed Description

[ImageHelper](#) (internal class, not intended for user level)

Helper for writing World images in DICOM. DICOM has a 'template' approach to image where MR [Image](#) Storage are distinct object from Enhanced MR [Image](#) Storage. For example the Pixel [Spacing](#) in one object is not at the same position (ie [Tag](#)) as in the other this class is the central (read: fragile) place where all the dispatching is done from a unified view of a world image (typically VTK or ITK point of view) down to the low level DICOM point of view.

#### Warning

: do not expect the API of this class to be maintained at any point, since as Modalities are added the API might have to be augmented or behavior changed to cope with new modalities.

### 27.146.2 Member Function Documentation

- 27.146.2.1 static [MediaStorage](#) [gdcm::ImageHelper::ComputeMediaStorageFromModality](#) ( const char \* *modality*, unsigned int *dimension* = 2, [PixelFormat](#) const & *pf* = [PixelFormat](#) (), [PhotometricInterpretation](#) const & *pi* = [PhotometricInterpretation](#) (), double *rescaleintercept* = 0, double *rescaleslope* = 1 ) [static]

Moved from [MediaStorage](#) here, since we need extra info stored in [PixelFormat](#) & [PhotometricInterpretation](#).

27.146.2.2 `static bool gdcm::ImageHelper::ComputeSpacingFromImagePositionPatient ( const std::vector< double > & imageposition, std::vector< double > & spacing ) [static]`

DO NOT USE.

27.146.2.3 `static std::vector<unsigned int> gdcm::ImageHelper::GetDimensionsValue ( const File & f ) [static]`

This function checks tags (0x0028, 0x0010) and (0x0028, 0x0011) for the rows and columns of the image in pixels (as opposed to actual distances). The output is {col , row}

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

27.146.2.4 `static bool gdcm::ImageHelper::GetDirectionCosinesFromDataSet ( DataSet const & ds, std::vector< double > & dircos ) [static]`

27.146.2.5 `static std::vector<double> gdcm::ImageHelper::GetDirectionCosinesValue ( File const & f ) [static]`

Get Direction Cosines (IOP) from/to a file Requires a file because mediastorage must be known

27.146.2.6 `static bool gdcm::ImageHelper::GetForcePixelSpacing ( ) [static]`

27.146.2.7 `static bool gdcm::ImageHelper::GetForceRescaleInterceptSlope ( ) [static]`

27.146.2.8 `static SmartPointer<LookupTable> gdcm::ImageHelper::GetLUT ( File const & f ) [static]`

27.146.2.9 `static std::vector<double> gdcm::ImageHelper::GetOriginValue ( File const & f ) [static]`

Set/Get Origin (IPP) from/to a file.

27.146.2.10 `static PhotometricInterpretation gdcm::ImageHelper::GetPhotometricInterpretationValue ( File const & f ) [static]`

27.146.2.11 `static PixelFormat gdcm::ImageHelper::GetPixelFormatValue ( const File & f ) [static]`

This function returns pixel information about an image from its dataset That includes samples per pixel and bit depth (in that order)

27.146.2.12 `static unsigned int gdcm::ImageHelper::GetPlanarConfigurationValue ( const File & f ) [static]`

27.146.2.13 `static const ByteValue* gdcm::ImageHelper::GetPointerFromElement ( Tag const & tag, File const & f ) [static]`

Moved from PixampReader to here. Generally used for photometric interpretation.

27.146.2.14 `static std::vector<double> gdcm::ImageHelper::GetRescaleInterceptSlopeValue ( File const & f ) [static]`

Set/Get shift/scale from/to a file

## Warning

this function reads/sets the Slope/Intercept in appropriate class storage, but also Grid Scaling in RT Dose Storage  
Can't take a dataset because the mediastorage of the file must be known

27.146.2.15 `static Tag gdcM::ImageHelper::GetSpacingTagFromMediaStorage ( MediaStorage const & ms ) [static],  
[protected]`

27.146.2.16 `static std::vector<double> gdcM::ImageHelper::GetSpacingValue ( File const & f ) [static]`

Set/Get [Spacing](#) from/to a [File](#).

27.146.2.17 `static Tag gdcM::ImageHelper::GetZSpacingTagFromMediaStorage ( MediaStorage const & ms ) [static],  
[protected]`

27.146.2.18 `static void gdcM::ImageHelper::SetDimensionsValue ( File & f, const Image & img ) [static]`

27.146.2.19 `static void gdcM::ImageHelper::SetDirectionCosinesValue ( DataSet & ds, const std::vector< double > & dircos )  
[static]`

Set Direction Cosines (IOP) from/to a file When [IOD](#) does not defines what is IOP (eg. typically Secondary Capture [Image](#) Storage) this call will simply remove the IOP attribute. Else in case of MR/CT image storage, this call will properly lookup the correct attribute to store the IOP.

27.146.2.20 `static void gdcM::ImageHelper::SetForcePixelSpacing ( bool ) [static]`

GDCM 1.x compatibility issue: When using ReWrite an MR [Image](#) Storage would be rewritten as Secondary Capture [Object](#) while still having a Pixel [Spacing](#) tag (0028,0030). If you have deal with those files, use this very special flag to handle them Unless explicitly set elsewhere by the standard, it will use value from 0028,0030 / 0018,0088 for the Pixel [Spacing](#) of the [Image](#)

27.146.2.21 `static void gdcM::ImageHelper::SetForceRescaleInterceptSlope ( bool ) [static]`

GDCM 1.x compatibility issue: when using ReWrite an MR [Image](#) Storage would be rewritten with a Rescale Slope/↵ Intercept while the standard would prohibit this (Philips Medical [System](#) is still doing that) Unless explicitly set elsewhere by the standard, it will use value from 0028,1052 / 0028,1053 for the Rescale Slope & Rescale Intercept values

27.146.2.22 `static void gdcM::ImageHelper::SetOriginValue ( DataSet & ds, const Image & img ) [static]`

27.146.2.23 `static void gdcM::ImageHelper::SetRescaleInterceptSlopeValue ( File & f, const Image & img ) [static]`

27.146.2.24 `static void gdcM::ImageHelper::SetSpacingValue ( DataSet & ds, const std::vector< double > & spacing )  
[static]`

The documentation for this class was generated from the following file:


- [gdcMImageHelper.h](#)

## 27.147 gdcm::ImageReader Class Reference


[ImageReader](#).

```
#include <gdcmImageReader.h>
```

Inheritance diagram for gdcm::ImageReader:

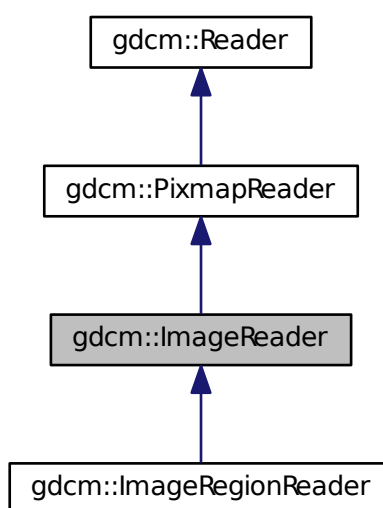

Collaboration diagram for `gdcm::ImageReader`:

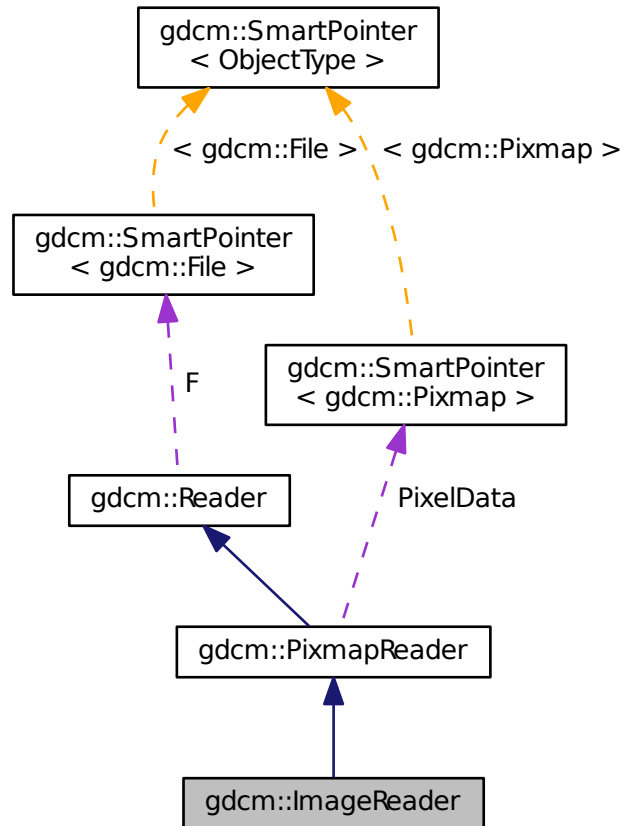

## Public Member Functions

- `ImageReader ()`
- `virtual ~ImageReader ()`
- `const Image & GetImage () const`  
*Return the read image.*
- `Image & GetImage ()`
- `virtual bool Read ()`

## Protected Member Functions

- `bool ReadACRNEMAIImage ()`
- `bool ReadImage (MediaStorage const &ms)`



## Additional Inherited Members

### 27.147.1 Detailed Description

[ImageReader](#).

#### Note

its role is to convert the DICOM [DataSet](#) into a [Image](#) representation [Image](#) is different from [Pixmap](#) has it has a position and a direction in Space.

#### See also

[Image](#)

#### Examples:

[CheckBigEndianBug.cxx](#), [CompressImage.cxx](#), [ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetJPEGSamplePrecision.cxx](#), [HelloVizWorld.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

### 27.147.2 Constructor & Destructor Documentation

27.147.2.1 `gdcm::ImageReader::ImageReader ( )`

27.147.2.2 `virtual gdcm::ImageReader::~~ImageReader ( ) [virtual]`

### 27.147.3 Member Function Documentation

27.147.3.1 `const Image& gdcm::ImageReader::GetImage ( ) const`

Return the read image.

#### Examples:

[CompressImage.cxx](#), [ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetJPEGSamplePrecision.cxx](#), [HelloVizWorld.cxx](#), [MergeTwoFiles.cxx](#), [PatchFile.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

27.147.3.2 `Image& gdcm::ImageReader::GetImage ( )`

27.147.3.3 `virtual bool gdcm::ImageReader::Read ( ) [virtual]`

Read the DICOM image. There are two reason for failure:

1. The input filename is not DICOM
2. The input DICOM file does not contains an [Image](#).

Reimplemented from [gdcm::PixmapReader](#).

Reimplemented in [gdcm::ImageRegionReader](#).

**Examples:**

[CheckBigEndianBug.cxx](#), [CompressImage.cxx](#), [ConvertToQImage.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetJPEGSamplePrecision.cxx](#), [HelloVizWorld.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [ReadMultiTimesException.cxx](#), and [threadgdcm.cxx](#).

**27.147.3.4** `bool gdcm::ImageReader::ReadACRNEMAIImage ( ) [protected],[virtual]`

Reimplemented from [gdcm::PixmapReader](#).

**27.147.3.5** `bool gdcm::ImageReader::ReadImage ( MediaStorage const & ms ) [protected],[virtual]`

Reimplemented from [gdcm::PixmapReader](#).

The documentation for this class was generated from the following file:


- [gdcmImageReader.h](#)

## 27.148 gdcm::ImageRegionReader Class Reference


[ImageRegionReader](#).

```
#include <gdcmImageRegionReader.h>
```

Inheritance diagram for `gdcm::ImageRegionReader`:

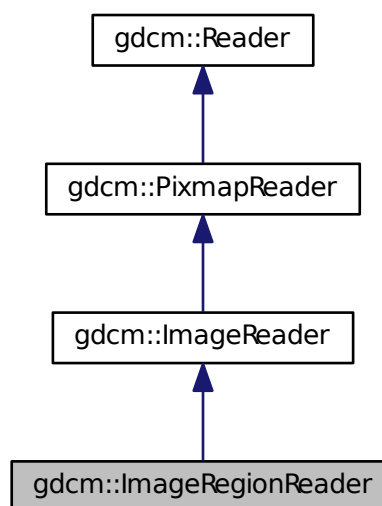

Collaboration diagram for gdcm::ImageRegionReader:

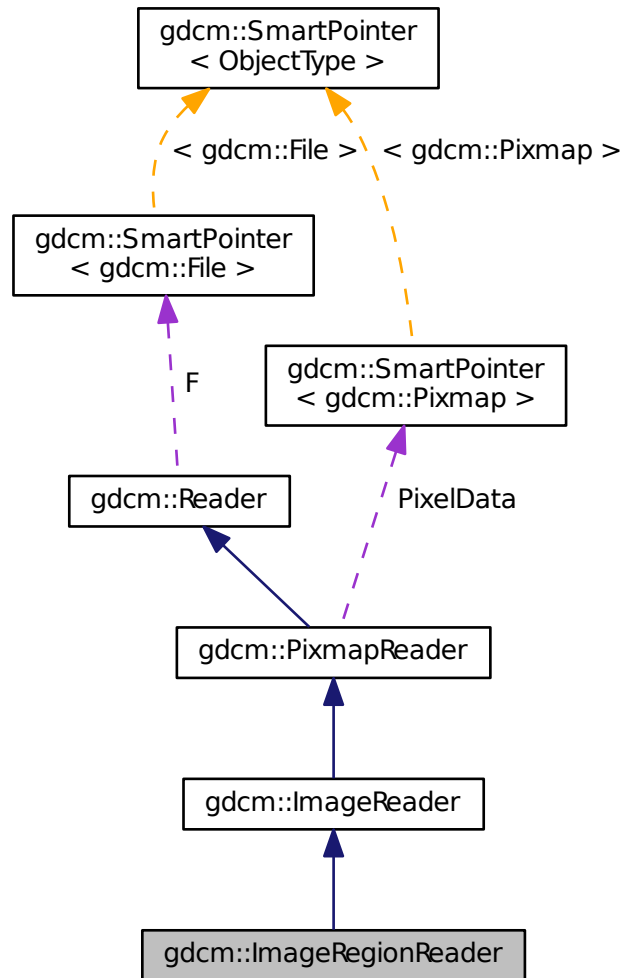

## Public Member Functions

- [ImageRegionReader](#) ()
- [~ImageRegionReader](#) ()
- [size\\_t ComputeBufferLength](#) () const
- [Region](#) const & [GetRegion](#) () const
- [bool ReadInformation](#) ()
- [bool ReadIntoBuffer](#) (char \*inreadbuffer, size\_t buflen)
- [void SetRegion](#) ([Region](#) const &region)

*Set/Get [Region](#) to be read.*

## Protected Member Functions

- bool [Read](#) ()

*To prevent user from calling super class [Read\(\)](#) function.*

## Additional Inherited Members

### 27.148.1 Detailed Description

[ImageRegionReader](#).

See also

[ImageReader](#)

### 27.148.2 Constructor & Destructor Documentation

27.148.2.1 `gdcm::ImageRegionReader::ImageRegionReader ( )`

27.148.2.2 `gdcm::ImageRegionReader::~~ImageRegionReader ( )`

### 27.148.3 Member Function Documentation

27.148.3.1 `size_t gdcm::ImageRegionReader::ComputeBufferLength ( ) const`

Explicit call which will compute the minimal buffer length that can hold the whole uncompressed image as defined by [Region](#) region.

Returns

0 upon error

27.148.3.2 `Region const& gdcm::ImageRegionReader::GetRegion ( ) const`

27.148.3.3 `bool gdcm::ImageRegionReader::Read ( ) [protected],[virtual]`

To prevent user from calling super class [Read\(\)](#) function.

Reimplemented from [gdcm::ImageReader](#).

27.148.3.4 `bool gdcm::ImageRegionReader::ReadInformation ( )`

Read meta information (not Pixel Data) from the DICOM file.

Returns

false upon error

27.148.3.5 bool gdcm::ImageRegionReader::ReadIntoBuffer ( char \* *inreadbuffer*, size\_t *buflen* )

Read into buffer:

Returns

false upon error

27.148.3.6 void gdcm::ImageRegionReader::SetRegion ( Region const & *region* )

Set/Get [Region](#) to be read.

The documentation for this class was generated from the following file:


- [gdcmImageRegionReader.h](#)

## 27.149 gdcm::ImageToImageFilter Class Reference


[ImageToImageFilter](#) class Super class for all filter taking an image and producing an output image.

```
#include <gdcmImageToImageFilter.h>
```

Inheritance diagram for gdcm::ImageToImageFilter:

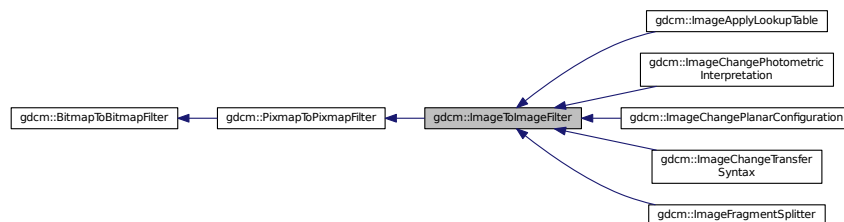

Collaboration diagram for `gdcm::ImageToImageFilter`:

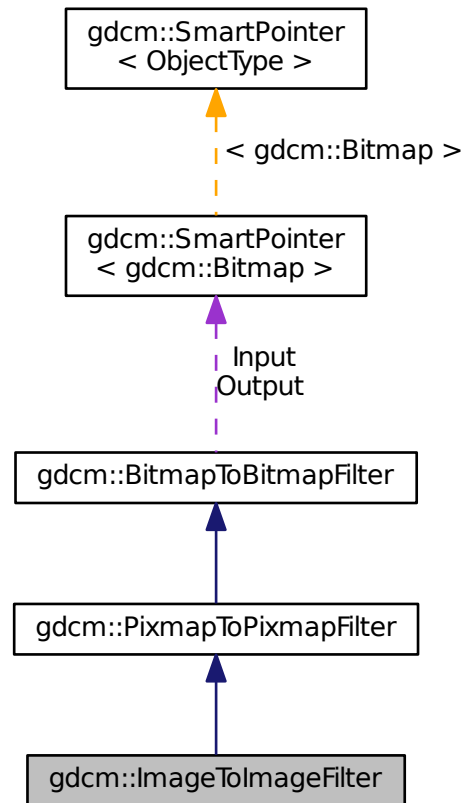

## Public Member Functions

- [ImageToImageFilter](#) ()
- [~ImageToImageFilter](#) ()
- [Image](#) & [GetInput](#) ()
- [const Image](#) & [GetOutput](#) () const  
*Get Output image.*

## Additional Inherited Members

### 27.149.1 Detailed Description

[ImageToImageFilter](#) class Super class for all filter taking an image and producing an output image.

### 27.149.2 Constructor & Destructor Documentation

27.149.2.1 `gdcm::ImageToImageFilter::ImageToImageFilter ( )`

27.149.2.2 `gdcm::ImageToImageFilter::~~ImageToImageFilter ( )` `[inline]`

### 27.149.3 Member Function Documentation

27.149.3.1 `Image& gdcm::ImageToImageFilter::GetInput ( )`

27.149.3.2 `const Image& gdcm::ImageToImageFilter::GetOutput ( ) const`

Get Output image.

Examples:

[CompressImage.cxx](#).

The documentation for this class was generated from the following file:


- [gdcmImageToImageFilter.h](#)

## 27.150 gdcm::ImageWriter Class Reference


[ImageWriter](#).

```
#include <gdcmImageWriter.h>
```

Inheritance diagram for `gdcm::ImageWriter`:

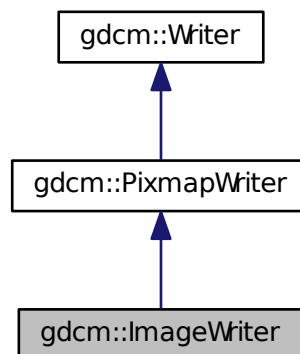

Collaboration diagram for `gdcm::ImageWriter`:

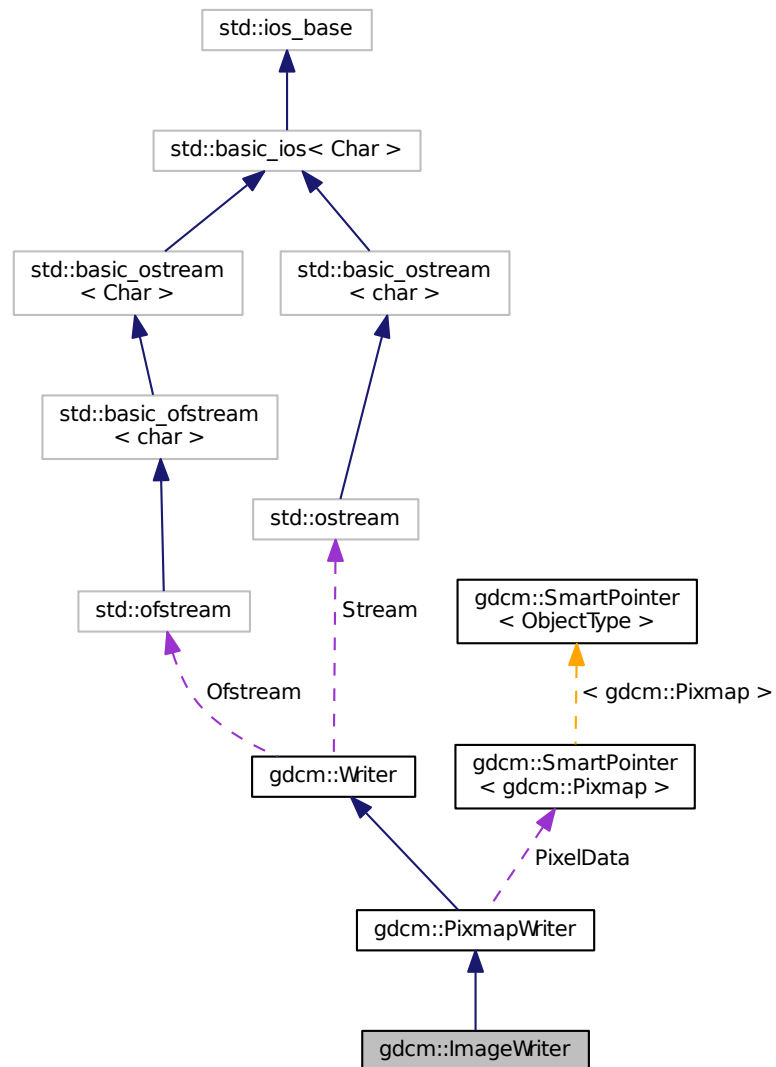

## Public Member Functions

- `ImageWriter ()`
- `~ImageWriter ()`
- `const Image & GetImage () const`
- `Image & GetImage ()`
- `bool Write ()`

*Write.*



## Additional Inherited Members

### 27.150.1 Detailed Description

[ImageWriter](#).

Examples:

[CompressImage.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [GenFakeImage.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), [iU22tomultisc.cxx](#), and [MergeTwoFiles.cxx](#).

### 27.150.2 Constructor & Destructor Documentation

27.150.2.1 `gdcm::ImageWriter::ImageWriter ( )`

27.150.2.2 `gdcm::ImageWriter::~~ImageWriter ( )`

### 27.150.3 Member Function Documentation

27.150.3.1 `const Image& gdcm::ImageWriter::GetImage ( ) const` `[inline],[virtual]`

Set/Get [Image](#) to be written It will overwrite anything [Image](#) infos found in [DataSet](#) (see parent class to see how to pass dataset)

Reimplemented from [gdcm::PixmapWriter](#).

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), and [iU22tomultisc.cxx](#).

27.150.3.2 `Image& gdcm::ImageWriter::GetImage ( )` `[inline],[virtual]`

Reimplemented from [gdcm::PixmapWriter](#).

27.150.3.3 `bool gdcm::ImageWriter::Write ( )` `[virtual]`

Write.

Reimplemented from [gdcm::Writer](#).

Examples:

[CompressImage.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [GenFakeImage.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), [iU22tomultisc.cxx](#), and [MergeTwoFiles.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmImageWriter.h](#)

## 27.151 gdcm::network::ImplementationClassUIDSub Class Reference

[ImplementationClassUIDSub](#) PS 3.7 [Table D.3-1](#) IMPLEMENTATION CLASS UID SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

```
#include <gdcmImplementationClassUIDSub.h>
```

### Public Member Functions

- [ImplementationClassUIDSub](#) ()
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.151.1 Detailed Description

[ImplementationClassUIDSub](#) PS 3.7 [Table D.3-1](#) IMPLEMENTATION CLASS UID SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

### 27.151.2 Constructor & Destructor Documentation

27.151.2.1 `gdcm::network::ImplementationClassUIDSub::ImplementationClassUIDSub ( )`

### 27.151.3 Member Function Documentation

27.151.3.1 `void gdcm::network::ImplementationClassUIDSub::Print ( std::ostream & os ) const`

27.151.3.2 `std::istream& gdcm::network::ImplementationClassUIDSub::Read ( std::istream & is )`

27.151.3.3 `size_t gdcm::network::ImplementationClassUIDSub::Size ( ) const`

27.151.3.4 `const std::ostream& gdcm::network::ImplementationClassUIDSub::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:

- [gdcmImplementationClassUIDSub.h](#)

## 27.152 gdcm::network::ImplementationUIDSub Class Reference

[ImplementationUIDSub](#) [Table D.3-2](#) IMPLEMENTATION UID SUB-ITEM FIELDS (A-ASSOCIATE-AC)

```
#include <gdcmImplementationUIDSub.h>
```

### Public Member Functions

- [ImplementationUIDSub](#) ()
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.152.1 Detailed Description

[ImplementationUIDSub Table](#) D.3-2 IMPLEMENTATION UID SUB-ITEM FIELDS (A-ASSOCIATE-AC)

### 27.152.2 Constructor & Destructor Documentation

27.152.2.1 `gdcm::network::ImplementationUIDSub::ImplementationUIDSub ( )`

### 27.152.3 Member Function Documentation

27.152.3.1 `const std::ostream& gdcm::network::ImplementationUIDSub::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:

- [gdcmImplementationUIDSub.h](#)

## 27.153 gdcm::network::ImplementationVersionNameSub Class Reference

[ImplementationVersionNameSub Table](#) D.3-3 IMPLEMENTATION VERSION NAME SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

```
#include <gdcmImplementationVersionNameSub.h>
```

### Public Member Functions

- [ImplementationVersionNameSub](#) ()
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.153.1 Detailed Description

[ImplementationVersionNameSub Table](#) D.3-3 IMPLEMENTATION VERSION NAME SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

### 27.153.2 Constructor & Destructor Documentation

27.153.2.1 `gdcm::network::ImplementationVersionNameSub::ImplementationVersionNameSub ( )`

### 27.153.3 Member Function Documentation

27.153.3.1 `void gdcm::network::ImplementationVersionNameSub::Print ( std::ostream & os ) const`

27.153.3.2 `std::istream& gdcm::network::ImplementationVersionNameSub::Read ( std::istream & is )`

27.153.3.3 `size_t gdcm::network::ImplementationVersionNameSub::Size ( ) const`

27.153.3.4 `const std::ostream& gdcmm::network::ImplementationVersionNameSub::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:


- [gdcmmImplementationVersionNameSub.h](#)

## 27.154 gdcmm::ImplicitDataElement Class Reference


Class to represent an *Implicit VR Data Element*.

```
#include <gdcmmImplicitDataElement.h>
```

Inheritance diagram for gdcmm::ImplicitDataElement:

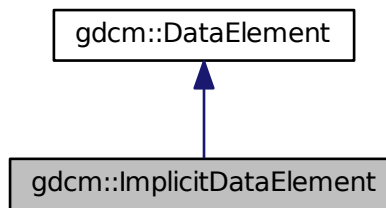

Collaboration diagram for gdcmm::ImplicitDataElement:

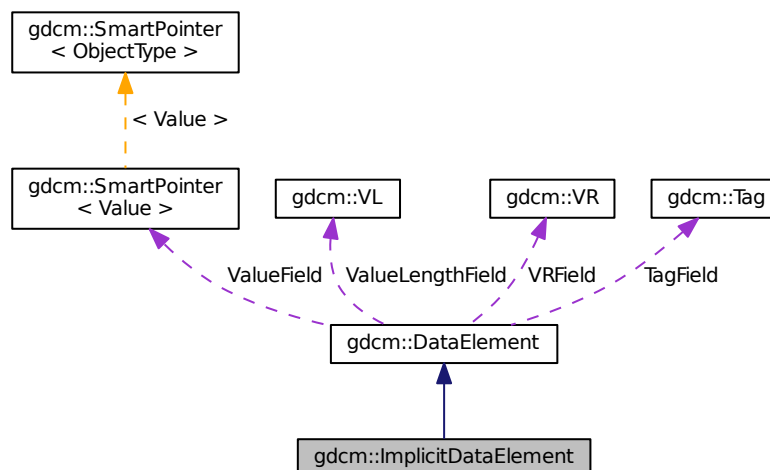

## Public Member Functions

- [VL GetLength](#) () const
- template<typename TSwap > std::istream & [Read](#) (std::istream &is)
- template<typename TSwap > std::istream & [ReadPreValue](#) (std::istream &is)
- template<typename TSwap > std::istream & [ReadValue](#) (std::istream &is, bool readvalues=true)
- template<typename TSwap > std::istream & [ReadValueWithLength](#) (std::istream &is, [VL](#) &length, bool readvalues=true)
- template<typename TSwap > std::istream & [ReadWithLength](#) (std::istream &is, [VL](#) &length, bool readvalues=true)
- template<typename TSwap > const std::ostream & [Write](#) (std::ostream &os) const

## Additional Inherited Members

### 27.154.1 Detailed Description

Class to represent an *Implicit VR Data Element*.

Note

bla

Examples:

[ReadExplicitLengthSQIVR.cxx](#).

### 27.154.2 Member Function Documentation

27.154.2.1 [VL gdcm::ImplicitDataElement::GetLength \( \) const](#)

27.154.2.2 [template<typename TSwap > std::istream& gdcm::ImplicitDataElement::Read \( std::istream & is \)](#)

27.154.2.3 [template<typename TSwap > std::istream& gdcm::ImplicitDataElement::ReadPreValue \( std::istream & is \)](#)

27.154.2.4 [template<typename TSwap > std::istream& gdcm::ImplicitDataElement::ReadValue \( std::istream & is, bool readvalues = true \)](#)

27.154.2.5 [template<typename TSwap > std::istream& gdcm::ImplicitDataElement::ReadValueWithLength \( std::istream & is, VL & length, bool readvalues = true \)](#)

27.154.2.6 [template<typename TSwap > std::istream& gdcm::ImplicitDataElement::ReadWithLength \( std::istream & is, VL & length, bool readvalues = true \)](#)

27.154.2.7 [template<typename TSwap > const std::ostream& gdcm::ImplicitDataElement::Write \( std::ostream & os \) const](#)


The documentation for this class was generated from the following file:

- [gdcmImplicitDataElement.h](#)


## 27.155 gdcm::InitializeEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for gdcm::InitializeEvent:

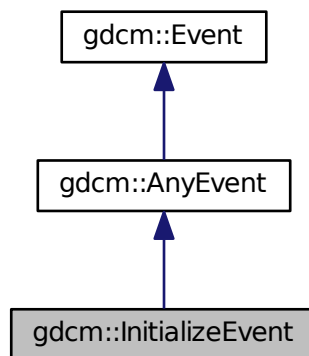

Collaboration diagram for gdcm::InitializeEvent:

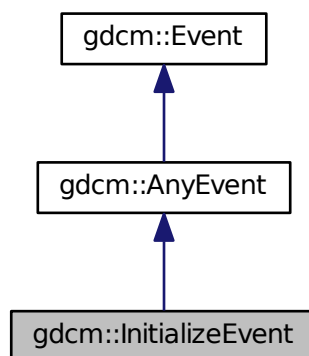

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.156 gdcm::IOD Class Reference

Class for representing a [IOD](#).

```
#include <gdcmIOD.h>
```

### Public Types

- typedef std::vector< [IODEntry](#) > [MapIODEntry](#)
- typedef MapIODEntry::size\_type [SizeType](#)

### Public Member Functions

- [IOD](#) ()
- void [AddIODEntry](#) (const [IODEntry](#) &iode)
- void [Clear](#) ()
- const [IODEntry](#) & [GetIODEntry](#) ([SizeType](#) idx) const
- [SizeType](#) [GetNumberOfIODs](#) () const
- [Type](#) [GetTypeFromTag](#) (const [Defs](#) &defs, const [Tag](#) &tag) const

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [IOD](#) &\_val)

### 27.156.1 Detailed Description

Class for representing a [IOD](#).

#### Note

bla

#### See also

[Dict](#)

#### Examples:

[TraverseModules.cxx](#).

### 27.156.2 Member Typedef Documentation

27.156.2.1 typedef std::vector<[IODEntry](#)> [gdcm::IOD::MapIODEntry](#)

27.156.2.2 typedef MapIODEntry::size\_type [gdcm::IOD::SizeType](#)

### 27.156.3 Constructor & Destructor Documentation

27.156.3.1 [gdcm::IOD::IOD](#) ( ) `[inline]`

### 27.156.4 Member Function Documentation

27.156.4.1 `void gdcmm::IOD::AddIODEntry ( const IODEntry & iode )` `[inline]`

27.156.4.2 `void gdcmm::IOD::Clear ( )` `[inline]`

27.156.4.3 `const IODEntry& gdcmm::IOD::GetIODEntry ( SizeType idx ) const` `[inline]`

Examples:

[TraverseModules.cxx](#).

27.156.4.4 `SizeType gdcmm::IOD::GetNumberOfIODs ( ) const` `[inline]`

Examples:

[TraverseModules.cxx](#).

27.156.4.5 `Type gdcmm::IOD::GetTypeFromTag ( const Defs & defs, const Tag & tag ) const`

### 27.156.5 Friends And Related Function Documentation

27.156.5.1 `std::ostream& operator<< ( std::ostream & _os, const IOD & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmmIOD.h](#)

## 27.157 gdcmm::IODEntry Class Reference

Class for representing a [IODEntry](#).

```
#include <gdcmmIODEntry.h>
```

### Public Member Functions

- [IODEntry](#) (const char \**name*="", const char \**ref*="", const char \**usag*="")
- const char \* [GetIE](#) () const
- const char \* [GetName](#) () const
- const char \* [GetRef](#) () const
- const char \* [GetUsage](#) () const
- [Usage::UsageType](#) [GetUsageType](#) () const
- void [SetIE](#) (const char \**ie*)
- void [SetName](#) (const char \**name*)
- void [SetRef](#) (const char \**ref*)
- void [SetUsage](#) (const char \**usag*)



## Friends

- `std::ostream & operator<< (std::ostream &_os, const IODEntry &_val)`

### 27.157.1 Detailed Description

Class for representing a [IODEntry](#).

#### Note

A.1.3 [IOD Module Table](#) and Functional Group [Macro Table](#) This Section of each [IOD](#) defines in a tabular form the [Modules](#) comprising the [IOD](#). The following information must be specified for each [Module](#) in the table:

- The name of the [Module](#) or Functional Group
- A reference to the Section in Annex C which defines the [Module](#) or Functional Group
- The usage of the [Module](#) or Functional Group; whether it is:
  - Mandatory (see A.1.3.1) , abbreviated M
  - Conditional (see A.1.3.2) , abbreviated C
  - User Option (see A.1.3.3) , abbreviated U
- The [Modules](#) referenced are defined in Annex C. A.1.3.1 MANDATORY MODULES For each [IOD](#), Mandatory [Modules](#) shall be supported per the definitions, semantics and requirements defined in Annex C. PS 3.3 - 2008 Page 96
- Standard - A.1.3.2 CONDITIONAL MODULES Conditional [Modules](#) are Mandatory [Modules](#) if specific conditions are met. If the specified conditions are not met, this [Module](#) shall not be supported; that is, no information defined in that [Module](#) shall be sent. A.1.3.3 USER OPTION MODULES User Option [Modules](#) may or may not be supported. If an optional [Module](#) is supported, the [Attribute](#) Types specified in the [Modules](#) in Annex C shall be supported.

#### See also

[DictEntry](#)

#### Examples:

[TraverseModules.cxx](#).

### 27.157.2 Constructor & Destructor Documentation

27.157.2.1 `gdcm::IODEntry::IODEntry ( const char * name = " ", const char * ref = " ", const char * usag = " " ) [inline]`

### 27.157.3 Member Function Documentation

27.157.3.1 `const char* gdcm::IODEntry::GetIE ( ) const [inline]`

27.157.3.2 `const char* gdcm::IODEntry::GetName ( ) const [inline]`

27.157.3.3 `const char* gdcm::IODEntry::GetRef ( ) const [inline]`

#### Examples:

[TraverseModules.cxx](#).

27.157.3.4 `const char* gdcm::IODEntry::GetUsage ( ) const` `[inline]`

27.157.3.5 `Usage::UsageType gdcm::IODEntry::GetUsageType ( ) const`

27.157.3.6 `void gdcm::IODEntry::SetIE ( const char * ie )` `[inline]`

27.157.3.7 `void gdcm::IODEntry::SetName ( const char * name )` `[inline]`

27.157.3.8 `void gdcm::IODEntry::SetRef ( const char * ref )` `[inline]`

27.157.3.9 `void gdcm::IODEntry::SetUsage ( const char * usag )` `[inline]`

## 27.157.4 Friends And Related Function Documentation

27.157.4.1 `std::ostream& operator<< ( std::ostream & _os, const IODEntry & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmIODEntry.h](#)

## 27.158 gdcm::IODs Class Reference

Class for representing a [IODs](#).

```
#include <gdcmIODs.h>
```

### Public Types

- typedef std::map< [IODName](#), [IOD](#) > [IODMapType](#)
- typedef IODMapType::const\_iterator [IODMapTypeConstIterator](#)
- typedef std::string [IODName](#)

### Public Member Functions

- [IODs](#) ()
- void [AddIOD](#) (const char \*name, const [IOD](#) &module)
- [IODMapTypeConstIterator](#) [Begin](#) () const
- void [Clear](#) ()
- [IODMapTypeConstIterator](#) [End](#) () const
- const [IOD](#) & [GetIOD](#) (const char \*name) const

### Friends

- std::ostream & [operator<<](#) (std::ostream & *\_os*, const [IODs](#) & *\_val*)

### 27.158.1 Detailed Description

Class for representing a [IODs](#).

Note

bla

See also

[IOD](#)

Examples:

[TraverseModules.cxx](#).

### 27.158.2 Member Typedef Documentation

27.158.2.1 `typedef std::map<IODName, IOD> gdcm::IODs::IODMapType`

27.158.2.2 `typedef IODMapType::const_iterator gdcm::IODs::IODMapTypeConstIterator`

27.158.2.3 `typedef std::string gdcm::IODs::IODName`

### 27.158.3 Constructor & Destructor Documentation

27.158.3.1 `gdcm::IODs::IODs ( )` `[inline]`

### 27.158.4 Member Function Documentation

27.158.4.1 `void gdcm::IODs::AddIOD ( const char * name, const IOD & module )` `[inline]`

27.158.4.2 `IODMapTypeConstIterator gdcm::IODs::Begin ( ) const` `[inline]`

Examples:

[TraverseModules.cxx](#).

27.158.4.3 `void gdcm::IODs::Clear ( )` `[inline]`

27.158.4.4 `IODMapTypeConstIterator gdcm::IODs::End ( ) const` `[inline]`

Examples:

[TraverseModules.cxx](#).

27.158.4.5 `const IOD& gdcm::IODs::GetIOD ( const char * name ) const` `[inline]`

### 27.158.5 Friends And Related Function Documentation

27.158.5.1 `std::ostream& operator<< ( std::ostream &_os, const IODs &_val )` [*friend*]

The documentation for this class was generated from the following file:


- [gdcmIODs.h](#)

## 27.159 gdcm::IPPSorter Class Reference


**IPPSorter** Implement a simple **Image** Position (**Patient**) sorter, along the **Image Orientation** (**Patient**) direction. This algorithm does NOT support duplicate and will FAIL in case of duplicate IPP.

```
#include <gdcmIPPSorter.h>
```

Inheritance diagram for `gdcm::IPPSorter`:

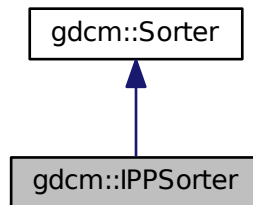

Collaboration diagram for `gdcm::IPPSorter`:

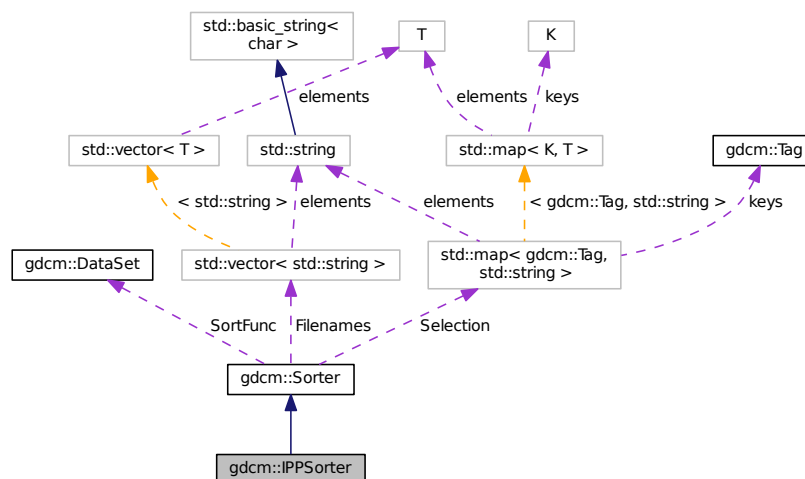

## Public Member Functions

- [IPPSorter](#) ()
- double [GetDirectionCosinesTolerance](#) () const
- double [GetZSpacing](#) () const
- double [GetZSpacingTolerance](#) () const
- void [SetComputeZSpacing](#) (bool b)
- void [SetDirectionCosinesTolerance](#) (double tol)
- void [SetDropDuplicatePositions](#) (bool b)
- void [SetZSpacingTolerance](#) (double tol)
- virtual bool [Sort](#) (std::vector< std::string > const &filenames)

## Protected Attributes

- bool [ComputeZSpacing](#)
- double [DirCosTolerance](#)
- bool [DropDuplicatePositions](#)
- double [ZSpacing](#)
- double [ZTolerance](#)

## Additional Inherited Members

### 27.159.1 Detailed Description

[IPPSorter](#) Implement a simple [Image](#) Position ([Patient](#)) sorter, along the [Image Orientation](#) ([Patient](#)) direction. This algorithm does NOT support duplicate and will FAIL in case of duplicate IPP.

#### Warning

See special note for [SetZSpacingTolerance](#) when computing the [ZSpacing](#) from the IPP of each DICOM files (default tolerance for consistent spacing is: 1e-6mm)

For more information on [Spacing](#), and how it is defined in DICOM, advanced users may refers to:

[http://gdcm.sourceforge.net/wiki/index.php/Imager\\_Pixel\\_Spacing](http://gdcm.sourceforge.net/wiki/index.php/Imager_Pixel_Spacing)

**Bug** There are currently a couple of bugs in this implementation:

- Gantry Tilt is not considered

#### Examples:

[Compute3DSpacing.cxx](#), [gdcmorthoplanes.cxx](#), [reslicesphere.cxx](#), and [VolumeSorter.cxx](#).

### 27.159.2 Constructor & Destructor Documentation

27.159.2.1 [gdcm::IPPSorter::IPPSorter](#) ( )

### 27.159.3 Member Function Documentation

27.159.3.1 [double gdcm::IPPSorter::GetDirectionCosinesTolerance](#) ( ) const `[inline]`

27.159.3.2 `double gdcmm::IPPSorter::GetZSpacing ( ) const [inline]`

Read-only function to provide access to the computed value for the Z-Spacing The ComputeZSpacing must have been set to true before execution of sort algorithm. Call this function *after* calling [Sort\(\)](#); Z-Spacing will be 0 on 2 occasions:

- Sorting simply failed, potentially duplicate IPP => ZSpacing = 0
- ZSpacing could not be computed (Z-Spacing is not constant, or ZTolerance is too low)

Examples:

[Compute3DSpacing.cxx](#), [gdcmmorthoplanes.cxx](#), and [reslicesphere.cxx](#).

27.159.3.3 `double gdcmm::IPPSorter::GetZSpacingTolerance ( ) const [inline]`

27.159.3.4 `void gdcmm::IPPSorter::SetComputeZSpacing ( bool b ) [inline]`

Functions related to Z-Spacing computation Set to true when sort algorithm should also perform a regular Z-Spacing computation using the [Image](#) Position ([Patient](#)) Potential reason for failure:

1. ALL slices are taken into account, if one slice is missing then ZSpacing will be set to 0 since the spacing will not be found to be regular along the [Series](#)

Examples:

[Compute3DSpacing.cxx](#), [gdcmmorthoplanes.cxx](#), [reslicesphere.cxx](#), and [VolumeSorter.cxx](#).

27.159.3.5 `void gdcmm::IPPSorter::SetDirectionCosinesTolerance ( double tol ) [inline]`

Sometimes IOP along a series is slightly changing for example: "0.999081\0.0426953\0.00369272\0.0419025\0.955059\0.293439", "0.999081\0.0426953\0.00369275\0.0419025\0.955059\0.293439", "0.999081\0.0426952\0.00369272\0.0419025\0.955059\0.293439", We need an API to define the tolerance which is allowed. Internally the cross vector of each direction cosines is computed. The tolerance then defines the distance in between 1. to the dot product of those cross vectors. In a perfect world this dot product is of course 1.0 which implies a [DirectionCosines](#) tolerance of exactly 0.0 (default).

27.159.3.6 `void gdcmm::IPPSorter::SetDropDuplicatePositions ( bool b ) [inline]`

Makes the [IPPSorter](#) ignore multiple images located at the same position. Only the first occurrence will be kept. DropDuplicatePositions defaults to false.

27.159.3.7 `void gdcmm::IPPSorter::SetZSpacingTolerance ( double tol ) [inline]`

1. Another reason for failure is that that Z-Spacing is only slightly changing (eg 1e-3) along the series, a human can determine that this is ok and change the tolerance from its default value: 1e-6

Examples:

[Compute3DSpacing.cxx](#), [gdcmmorthoplanes.cxx](#), and [reslicesphere.cxx](#).

27.159.3.8 `virtual bool gdcm::IPPSorter::Sort ( std::vector< std::string > const & filenames )` `[virtual]`

Main entry point to the sorter. It will execute the filter, option should be set before running this function (SetZSpacing↔Tolerance, ...) Return value indicate if sorting could be achived. Warning this does *NOT* imply that spacing is consistent, it only means the file are sorted according to IPP You should check if ZSpacing is 0 or not to deduce if file are actually a 3D volume

Reimplemented from [gdcm::Sorter](#).

Examples:

[Compute3DSpacing.cxx](#), [gdcmorthoplanes.cxx](#), [reslicesphere.cxx](#), and [VolumeSorter.cxx](#).

## 27.159.4 Member Data Documentation

27.159.4.1 `bool gdcm::IPPSorter::ComputeZSpacing` `[protected]`

27.159.4.2 `double gdcm::IPPSorter::DirCosTolerance` `[protected]`

27.159.4.3 `bool gdcm::IPPSorter::DropDuplicatePositions` `[protected]`

27.159.4.4 `double gdcm::IPPSorter::ZSpacing` `[protected]`

27.159.4.5 `double gdcm::IPPSorter::ZTolerance` `[protected]`

The documentation for this class was generated from the following file:


- [gdcmIPPSorter.h](#)

## 27.160 gdcm::Item Class Reference


Class to represent an [Item](#) A component of the value of a Data [Element](#) that is of [Value](#) Representation Sequence of Items. An [Item](#) contains a Data Set . See PS 3.5 7.5.1 [Item](#) Encoding Rules Each [Item](#) of a Data [Element](#) of VR SQ shall be encoded as a DICOM Standart Data [Element](#) with a specific Data [Element](#) Tag of [Value](#) (FFFE,E000). The [Item](#) [Tag](#) is followed by a 4 byte [Item](#) Length field encoded in one of the following two ways Explicit/ Implicit.

```
#include <gdcmItem.h>
```

Inheritance diagram for gdcm::Item:

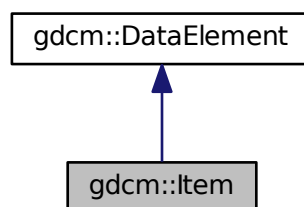

Collaboration diagram for `gdcm::Item`:

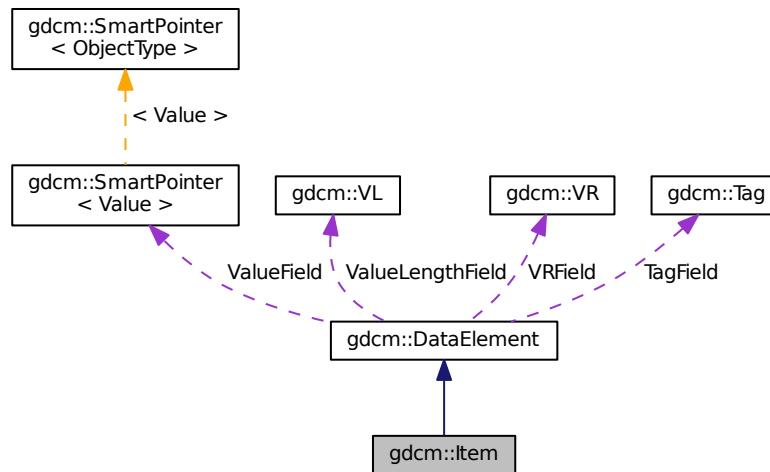

## Public Member Functions

- `Item ()`
- `Item (Item const &val)`
- `void Clear ()`
- `bool FindDataElement (const Tag &t) const`
- `const DataElement & GetDataElement (const Tag &t) const`
- `template<typename TDE > VL GetLength () const`
- `const DataSet & GetNestedDataSet () const`
- `DataSet & GetNestedDataSet ()`
- `void InsertDataElement (const DataElement &de)`
- `template<typename TDE , typename TSwap > std::istream & Read (std::istream &is)`
- `void SetNestedDataSet (const DataSet &nested)`
- `template<typename TDE , typename TSwap > const std::ostream & Write (std::ostream &os) const`

## Friends

- `std::ostream & operator<< (std::ostream &os, const Item &val)`

## Additional Inherited Members

### 27.160.1 Detailed Description

Class to represent an `Item`. A component of the value of a Data `Element` that is of `Value` Representation Sequence of Items. An `Item` contains a Data Set. See PS 3.5 7.5.1 `Item` Encoding Rules. Each `Item` of a Data `Element` of `VR` SQ



shall be encoded as a DICOM Standard Data [Element](#) with a specific Data [Element Tag](#) of [Value](#) (FFFE,E000). The [Item Tag](#) is followed by a 4 byte [Item Length](#) field encoded in one of the following two ways Explicit/ Implicit.

#### Note

ITEM: A component of the [Value](#) of a Data [Element](#) that is of [Value](#) Representation Sequence of Items. An [Item](#) contains a Data Set.

#### Examples:

[ChangeSequenceUltrasound.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllIVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), and [LargeVRDSExplicit.cxx](#).

## 27.160.2 Constructor & Destructor Documentation

27.160.2.1 `gdcm::Item::Item ( )` `[inline]`

27.160.2.2 `gdcm::Item::Item ( Item const & val )` `[inline]`

## 27.160.3 Member Function Documentation

27.160.3.1 `void gdcm::Item::Clear ( )` `[inline]`

References `gdcm::DataElement::Clear()`.

Referenced by `gdcm::SequenceOfItems::Read()`.

27.160.3.2 `bool gdcm::Item::FindDataElement ( const Tag & t ) const` `[inline]`

#### Examples:

[ReadAndDumpDICOMDIR.cxx](#).

27.160.3.3 `const DataElement& gdcm::Item::GetDataElement ( const Tag & t ) const` `[inline]`

#### Examples:

[ReadAndDumpDICOMDIR.cxx](#).

27.160.3.4 `template<typename TDE > VL gdcm::Item::GetLength ( ) const`

27.160.3.5 `const DataSet& gdcm::Item::GetNestedDataSet ( ) const` `[inline]`

#### Examples:

[ChangeSequenceUltrasound.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllIVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenSeqs.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), and [LargeVRDSExplicit.cxx](#).

Referenced by `gdcm::SequenceOfItems::Read()`.

27.160.3.6 `DataSet& gdcM::Item::GetNestedDataSet ( ) [inline]`

27.160.3.7 `void gdcM::Item::InsertDataElement ( const DataSet & de ) [inline]`

27.160.3.8 `template<typename TDE , typename TSwap > std::istream& gdcM::Item::Read ( std::istream & is ) [inline]`

References `gdcM::DataSet::Clear()`, `gdcMDebugMacro`, `gdcMErrorMacro`, `gdcMWarningMacro`, `gdcM::DataSet::IsEmpty()`, and `gdcM::SwapperDoOp::Swap()`.

Referenced by `gdcM::SequenceOfItems::Read()`.

27.160.3.9 `void gdcM::Item::SetNestedDataSet ( const DataSet & nested ) [inline]`

27.160.3.10 `template<typename TDE , typename TSwap > const std::ostream& gdcM::Item::Write ( std::ostream & os ) const [inline]`

References `gdcMWarningMacro`, `gdcM::VL::GetLength()`, `gdcM::VL::Write()`, and `gdcM::Tag::Write()`.

## 27.160.4 Friends And Related Function Documentation

27.160.4.1 `std::ostream& operator<< ( std::ostream & os, const Item & val ) [friend]`


The documentation for this class was generated from the following file:

- [gdcMItem.h](#)


## 27.161 gdcM::IterationEvent Class Reference

```
#include <gdcMEvent.h>
```

Inheritance diagram for `gdcM::IterationEvent`:

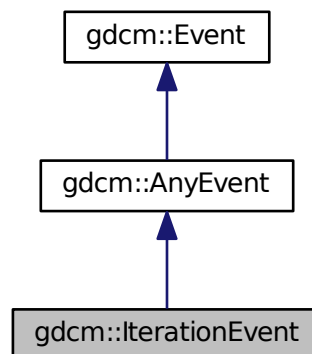

Collaboration diagram for gdcm::IterationEvent:

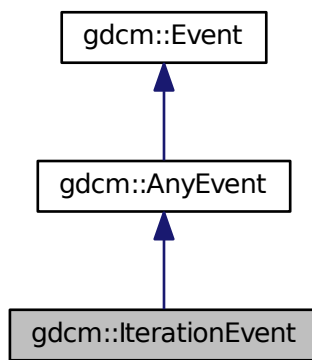

### Additional Inherited Members

The documentation for this class was generated from the following file:


- [gdcmEvent.h](#)

## 27.162 gdcm::JPEG12Codec Class Reference


Class to do JPEG 12bits (lossy & lossless)

```
#include <gdcmJPEG12Codec.h>
```

Inheritance diagram for `gdcm::JPEG12Codec`:

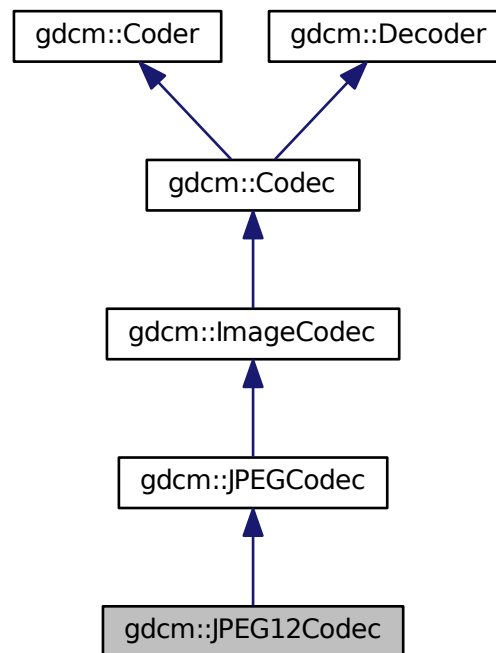

Collaboration diagram for `gdcm::JPEG12Codec`:

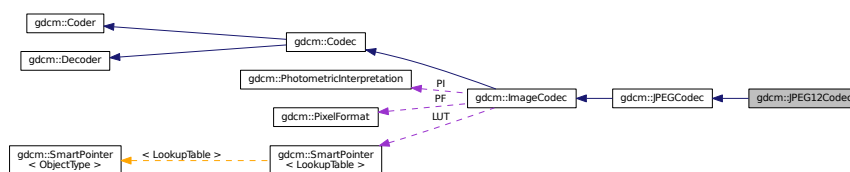

## Public Member Functions

- [JPEG12Codec](#) ()
- [~JPEG12Codec](#) ()
- [DecodeByStreams](#) (std::istream &is, std::ostream &os)
- [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- [InternalCode](#) (const char \*input, unsigned long len, std::ostream &os)

## Protected Member Functions

- virtual bool [EncodeBuffer](#) (std::ostream &os, const char \*data, size\_t datalen)
- bool [IsStateSuspension](#) () const

## Additional Inherited Members

### 27.162.1 Detailed Description

Class to do JPEG 12bits (lossy & lossless)

#### Note

internal class

### 27.162.2 Constructor & Destructor Documentation

27.162.2.1 `gdcm::JPEG12Codec::JPEG12Codec ( )`

27.162.2.2 `gdcm::JPEG12Codec::~~JPEG12Codec ( )`

### 27.162.3 Member Function Documentation

27.162.3.1 `bool gdcm::JPEG12Codec::DecodeByStreams ( std::istream & is, std::ostream & os )` [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.162.3.2 `virtual bool gdcm::JPEG12Codec::EncodeBuffer ( std::ostream & os, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.162.3.3 `bool gdcm::JPEG12Codec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.162.3.4 `bool gdcm::JPEG12Codec::InternalCode ( const char * input, unsigned long len, std::ostream & os )` [virtual]

Reimplemented from [gdcm::Coder](#).

27.162.3.5 `bool gdcm::JPEG12Codec::IsStateSuspension ( ) const` [protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

The documentation for this class was generated from the following file:


- [gdcmJPEG12Codec.h](#)

## 27.163 gdcm::JPEG16Codec Class Reference


Class to do JPEG 16bits (lossless)

```
#include <gdcmJPEG16Codec.h>
```

Inheritance diagram for gdcm::JPEG16Codec:

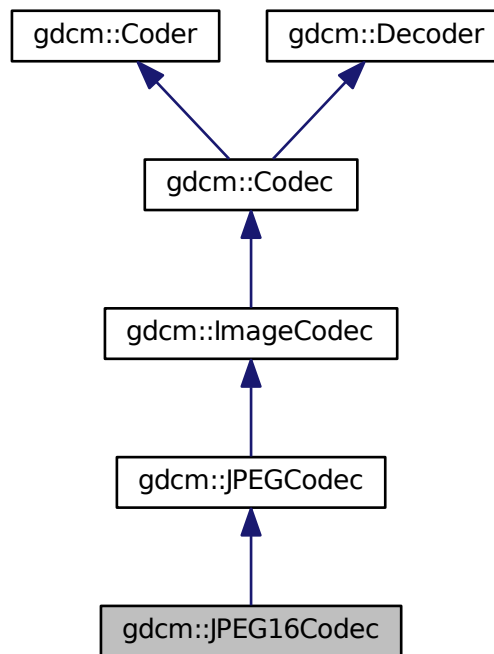

Collaboration diagram for gdcm::JPEG16Codec:

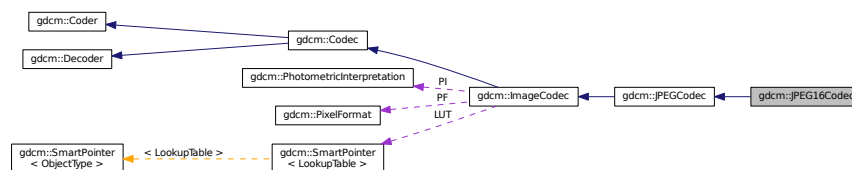

### Public Member Functions

- [JPEG16Codec](#) ()
- [~JPEG16Codec](#) ()
- [bool DecodeByStreams](#) (std::istream &is, std::ostream &os)

- bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- bool [InternalCode](#) (const char \*input, unsigned long len, std::ostream &os)

### Protected Member Functions

- virtual bool [EncodeBuffer](#) (std::ostream &os, const char \*data, size\_t datalen)
- bool [IsStateSuspension](#) () const

### Additional Inherited Members

#### 27.163.1 Detailed Description

Class to do JPEG 16bits (lossless)

##### Note

internal class

#### 27.163.2 Constructor & Destructor Documentation

27.163.2.1 [gdcm::JPEG16Codec::JPEG16Codec](#) ( )

27.163.2.2 [gdcm::JPEG16Codec::~~JPEG16Codec](#) ( )

#### 27.163.3 Member Function Documentation

27.163.3.1 [bool gdcm::JPEG16Codec::DecodeByStreams](#) ( std::istream & *is*, std::ostream & *os* ) [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.163.3.2 [virtual bool gdcm::JPEG16Codec::EncodeBuffer](#) ( std::ostream & *os*, const char \* *data*, size\_t *datalen* )  
[protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.163.3.3 [bool gdcm::JPEG16Codec::GetHeaderInfo](#) ( std::istream & *is*, [TransferSyntax](#) & *ts* ) [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.163.3.4 [bool gdcm::JPEG16Codec::InternalCode](#) ( const char \* *input*, unsigned long *len*, std::ostream & *os* ) [virtual]

Reimplemented from [gdcm::Coder](#).

27.163.3.5 [bool gdcm::JPEG16Codec::IsStateSuspension](#) ( ) const [protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

The documentation for this class was generated from the following file:


- [gdcmJPEG16Codec.h](#)

## 27.164 gdcm::JPEG2000Codec Class Reference


Class to do JPEG 2000.

```
#include <gdcmJPEG2000Codec.h>
```

Inheritance diagram for gdcm::JPEG2000Codec:

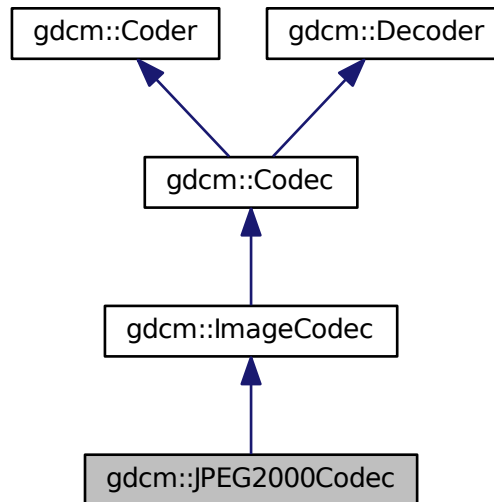

Collaboration diagram for gdcm::JPEG2000Codec:

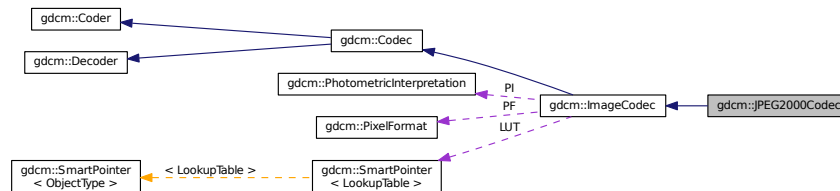

### Public Member Functions

- [JPEG2000Codec](#) ()
- [~JPEG2000Codec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &ts) const  
Return whether this coder support this transfer syntax (can code it)
- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const  
Return whether this decoder support this transfer syntax (can decode it)



- virtual [ImageCodec](#) \* [Clone](#) () const
- bool [Code](#) ([DataElement](#) const &in, [DataElement](#) &out)  
*Code.*
- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)  
*Decode.*
- virtual bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- double [GetQuality](#) (unsigned int idx=0) const
- double [GetRate](#) (unsigned int idx=0) const
- void [SetNumberOfResolutions](#) (unsigned int nres)
- void [SetQuality](#) (unsigned int idx, double q)
- void [SetRate](#) (unsigned int idx, double rate)
- void [SetReversible](#) (bool res)
- void [SetTileSize](#) (unsigned int tx, unsigned int ty)

### Protected Member Functions

- bool [AppendFrameEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [AppendRowEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [DecodeByStreams](#) (std::istream &is, std::ostream &os)
- bool [DecodeExtent](#) (char \*buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream &is)
- bool [IsFrameEncoder](#) ()
- bool [IsRowEncoder](#) ()
- bool [StartEncode](#) (std::ostream &)
- bool [StopEncode](#) (std::ostream &)

### Friends

- class [Bitmap](#)
- class [ImageRegionReader](#)

### Additional Inherited Members

#### 27.164.1 Detailed Description

Class to do JPEG 2000.

#### Note

the class will produce JPC (JPEG 2000 codestream), since some private implementor are using full jp2 file the decoder tolerate jp2 input this is an implementation of an [ImageCodec](#)

#### 27.164.2 Constructor & Destructor Documentation

27.164.2.1 [gdcm::JPEG2000Codec::JPEG2000Codec \( \)](#)

27.164.2.2 [gdcm::JPEG2000Codec::~~JPEG2000Codec \( \)](#)

### 27.164.3 Member Function Documentation

**27.164.3.1** `bool gdcm::JPEG2000Codec::AppendFrameEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.2** `bool gdcm::JPEG2000Codec::AppendRowEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.3** `bool gdcm::JPEG2000Codec::CanCode ( TransferSyntax const & ) const` [virtual]

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.4** `bool gdcm::JPEG2000Codec::CanDecode ( TransferSyntax const & ) const` [virtual]

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.5** `virtual ImageCodec* gdcm::JPEG2000Codec::Clone ( ) const` [virtual]

Implements [gdcm::ImageCodec](#).

**27.164.3.6** `bool gdcm::JPEG2000Codec::Code ( DataElement const & in_, DataElement & out_ )` [virtual]

Code.

Reimplemented from [gdcm::Coder](#).

**27.164.3.7** `bool gdcm::JPEG2000Codec::Decode ( DataElement const &, DataElement & )` [virtual]

Decode.

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.8** `bool gdcm::JPEG2000Codec::DecodeByStreams ( std::istream & is, std::ostream & os )` [protected],  
[virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.164.3.9** `bool gdcm::JPEG2000Codec::DecodeExtent ( char * buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin,  
unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream & is )` [protected]

27.164.3.10 `virtual bool gdcm::JPEG2000Codec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.164.3.11 `double gdcm::JPEG2000Codec::GetQuality ( unsigned int idx = 0 ) const`

27.164.3.12 `double gdcm::JPEG2000Codec::GetRate ( unsigned int idx = 0 ) const`

27.164.3.13 `bool gdcm::JPEG2000Codec::IsFrameEncoder ( )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.164.3.14 `bool gdcm::JPEG2000Codec::IsRowEncoder ( )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.164.3.15 `void gdcm::JPEG2000Codec::SetNumberOfResolutions ( unsigned int nres )`

27.164.3.16 `void gdcm::JPEG2000Codec::SetQuality ( unsigned int idx, double q )`

27.164.3.17 `void gdcm::JPEG2000Codec::SetRate ( unsigned int idx, double rate )`

27.164.3.18 `void gdcm::JPEG2000Codec::SetReversible ( bool res )`

27.164.3.19 `void gdcm::JPEG2000Codec::SetTileSize ( unsigned int tx, unsigned int ty )`

27.164.3.20 `bool gdcm::JPEG2000Codec::StartEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.164.3.21 `bool gdcm::JPEG2000Codec::StopEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

## 27.164.4 Friends And Related Function Documentation

27.164.4.1 `friend class Bitmap` [friend]

27.164.4.2 `friend class ImageRegionReader` [friend]

The documentation for this class was generated from the following file:


- [gdcmJPEG2000Codec.h](#)

## 27.165 gdcm::JPEG8Codec Class Reference


Class to do JPEG 8bits (lossy & lossless)

```
#include <gdcmJPEG8Codec.h>
```

Inheritance diagram for `gdcm::JPEG8Codec`:

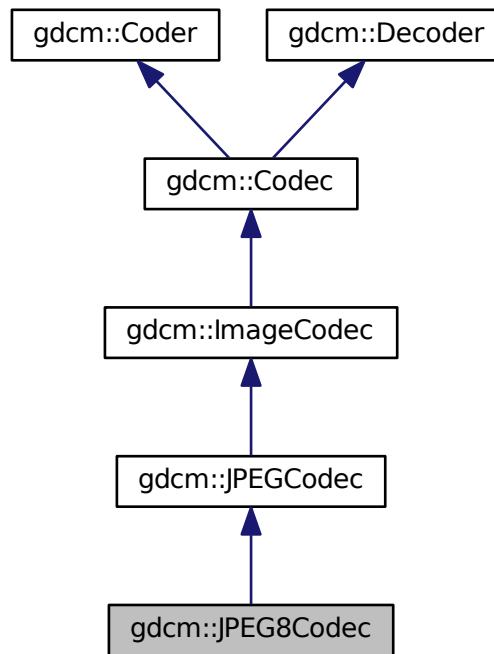

Collaboration diagram for `gdcm::JPEG8Codec`:

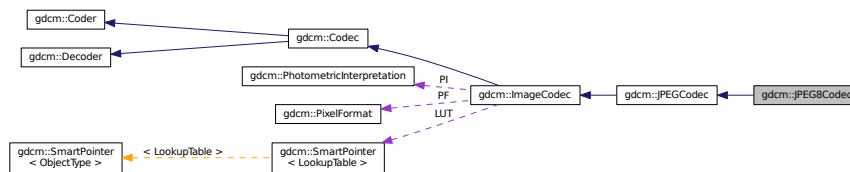

## Public Member Functions

- [JPEG8Codec](#) ()
- [~JPEG8Codec](#) ()
- [DecodeByStreams](#) (std::istream &is, std::ostream &os)
- [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- [InternalCode](#) (const char \*input, unsigned long len, std::ostream &os)

## Protected Member Functions

- virtual bool [EncodeBuffer](#) (std::ostream &os, const char \*data, size\_t datalen)
- bool [IsStateSuspension](#) () const

## Additional Inherited Members

### 27.165.1 Detailed Description

Class to do JPEG 8bits (lossy & lossless)

#### Note

internal class

### 27.165.2 Constructor & Destructor Documentation

27.165.2.1 `gdcm::JPEG8Codec::JPEG8Codec ( )`

27.165.2.2 `gdcm::JPEG8Codec::~~JPEG8Codec ( )`

### 27.165.3 Member Function Documentation

27.165.3.1 `bool gdcm::JPEG8Codec::DecodeByStreams ( std::istream & is, std::ostream & os )` [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.165.3.2 `virtual bool gdcm::JPEG8Codec::EncodeBuffer ( std::ostream & os, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.165.3.3 `bool gdcm::JPEG8Codec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

27.165.3.4 `bool gdcm::JPEG8Codec::InternalCode ( const char * input, unsigned long len, std::ostream & os )` [virtual]

Reimplemented from [gdcm::Coder](#).

27.165.3.5 `bool gdcm::JPEG8Codec::IsStateSuspension ( ) const` [protected], [virtual]

Reimplemented from [gdcm::JPEGCodec](#).

The documentation for this class was generated from the following file:


- [gdcmJPEG8Codec.h](#)

## 27.166 gdcm::JPEGCodec Class Reference


JPEG codec Class to do JPEG (8bits, 12bits, 16bits lossy & lossless). It redispatch in between the different codec implementation: [JPEG8Codec](#), [JPEG12Codec](#) & [JPEG16Codec](#). It also support inconsistency in between DICOM header and JPEG compressed stream [ImageCodec](#) implementation for the JPEG case.

```
#include <gdcmJPEGCodec.h>
```

Inheritance diagram for gdcm::JPEGCodec:

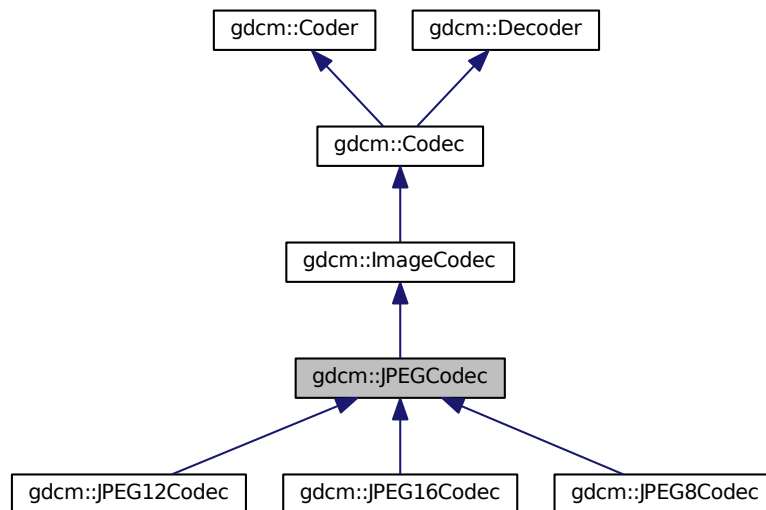

Collaboration diagram for gdcm::JPEGCodec:

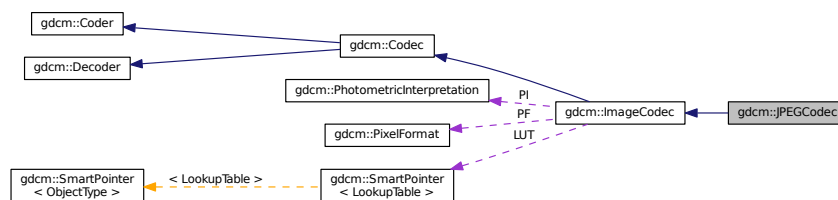

### Public Member Functions

- [JPEGCodec](#) ()
- [~JPEGCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &ts) const  
Return whether this coder support this transfer syntax (can code it)
- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const

*Return whether this decoder support this transfer syntax (can decode it)*

- virtual [ImageCodec](#) \* [Clone](#) () const
- bool [Code](#) ([DataElement](#) const &in, [DataElement](#) &out)

*Compress into JPEG.*

- void [ComputeOffsetTable](#) (bool b)

*Compute the offset table:*

- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)

*Decode.*

- virtual bool [EncodeBuffer](#) (std::ostream &out, const char \*inbuffer, size\_t inlen)
- virtual bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- bool [GetLossless](#) () const
- double [GetQuality](#) () const
- void [SetLossless](#) (bool l)
- void [SetPixelFormat](#) ([PixelFormat](#) const &pf)
- void [SetQuality](#) (double q)

## Protected Member Functions

- bool [AppendFrameEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [AppendRowEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [DecodeByStreams](#) (std::istream &is, std::ostream &os)
- bool [DecodeExtent](#) (char \*buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream &is)
- bool [IsFrameEncoder](#) ()
- bool [IsRowEncoder](#) ()
- virtual bool [IsStateSuspension](#) () const
- bool [IsValid](#) ([PhotometricInterpretation](#) const &pi)
- void [SetBitSample](#) (int bit)
- bool [StartEncode](#) (std::ostream &)
- bool [StopEncode](#) (std::ostream &)

## Protected Attributes

- int [BitSample](#)
- int [Quality](#)

## Friends

- class [ImageRegionReader](#)

## Additional Inherited Members

### 27.166.1 Detailed Description

JPEG codec Class to do JPEG (8bits, 12bits, 16bits lossy & lossless). It redispach in between the different codec implementation: [JPEG8Codec](#), [JPEG12Codec](#) & [JPEG16Codec](#) It also support inconsistency in between DICO↔M header and JPEG compressed stream [ImageCodec](#) implementation for the JPEG case.

## Note

Things you should know if you ever want to dive into DICOM/JPEG world (among other):

- [http://groups.google.com/group/comp.protocols.dicom/browse\\_thread/thread/625e46919f208](http://groups.google.com/group/comp.protocols.dicom/browse_thread/thread/625e46919f208)
- [http://groups.google.com/group/comp.protocols.dicom/browse\\_thread/thread/75fdfccc65a62](http://groups.google.com/group/comp.protocols.dicom/browse_thread/thread/75fdfccc65a62)
- [http://groups.google.com/group/comp.protocols.dicom/browse\\_thread/thread/2d525ef6a2f09](http://groups.google.com/group/comp.protocols.dicom/browse_thread/thread/2d525ef6a2f09)
- [http://groups.google.com/group/comp.protocols.dicom/browse\\_thread/thread/6b93af410f8c9](http://groups.google.com/group/comp.protocols.dicom/browse_thread/thread/6b93af410f8c9)

## Examples:

[GetJPEGSamplePrecision.cxx](#).

## 27.166.2 Constructor & Destructor Documentation

27.166.2.1 `gdcm::JPEGCodec::JPEGCodec ( )`

27.166.2.2 `gdcm::JPEGCodec::~~JPEGCodec ( )`

## 27.166.3 Member Function Documentation

27.166.3.1 `bool gdcm::JPEGCodec::AppendFrameEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.2 `bool gdcm::JPEGCodec::AppendRowEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.3 `bool gdcm::JPEGCodec::CanCode ( TransferSyntax const & ) const` [virtual]

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.4 `bool gdcm::JPEGCodec::CanDecode ( TransferSyntax const & ) const` [virtual]

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.5 `virtual ImageCodec* gdcm::JPEGCodec::Clone ( ) const` [virtual]

Implements [gdcm::ImageCodec](#).



27.166.3.6 `bool gdcm::JPEGCodec::Code ( DataElement const & in, DataElement & out )` [virtual]

Compress into JPEG.

Reimplemented from [gdcm::Coder](#).

27.166.3.7 `void gdcm::JPEGCodec::ComputeOffsetTable ( bool b )`

Compute the offset table:

27.166.3.8 `bool gdcm::JPEGCodec::Decode ( DataElement const & , DataElement & )` [virtual]

Decode.

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.9 `bool gdcm::JPEGCodec::DecodeByStreams ( std::istream & is, std::ostream & os )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.10 `bool gdcm::JPEGCodec::DecodeExtent ( char * buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream & is )` [protected]

27.166.3.11 `virtual bool gdcm::JPEGCodec::EncodeBuffer ( std::ostream & out, const char * inbuffer, size_t inlen )`  
[virtual]

Reimplemented in [gdcm::JPEG12Codec](#), [gdcm::JPEG16Codec](#), and [gdcm::JPEG8Codec](#).

27.166.3.12 `virtual bool gdcm::JPEGCodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented from [gdcm::ImageCodec](#).

Reimplemented in [gdcm::JPEG12Codec](#), [gdcm::JPEG16Codec](#), and [gdcm::JPEG8Codec](#).

Examples:

[GetJPEGSamplePrecision.cxx](#).

27.166.3.13 `bool gdcm::JPEGCodec::GetLossless ( )` const

27.166.3.14 `double gdcm::JPEGCodec::GetQuality ( )` const

27.166.3.15 `bool gdcm::JPEGCodec::IsFrameEncoder ( )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.16 `bool gdcm::JPEGCodec::IsRowEncoder ( )` [protected],[virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.166.3.17 `virtual bool gdcM::JPEGCodec::IsStateSuspension ( ) const` [protected],[virtual]

Reimplemented in [gdcM::JPEG12Codec](#), [gdcM::JPEG16Codec](#), and [gdcM::JPEG8Codec](#).

27.166.3.18 `bool gdcM::JPEGCodec::IsValid ( PhotometricInterpretation const & pi )` [protected],[virtual]

Reimplemented from [gdcM::ImageCodec](#).

27.166.3.19 `void gdcM::JPEGCodec::SetBitSample ( int bit )` [protected]

27.166.3.20 `void gdcM::JPEGCodec::SetLossless ( bool l )`

27.166.3.21 `void gdcM::JPEGCodec::SetPixelFormat ( PixelFormat const & pf )` [virtual]

Reimplemented from [gdcM::ImageCodec](#).

Examples:

[GetJPEGSamplePrecision.cxx](#).

27.166.3.22 `void gdcM::JPEGCodec::SetQuality ( double q )`

27.166.3.23 `bool gdcM::JPEGCodec::StartEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcM::ImageCodec](#).

27.166.3.24 `bool gdcM::JPEGCodec::StopEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcM::ImageCodec](#).

## 27.166.4 Friends And Related Function Documentation

27.166.4.1 `friend class ImageRegionReader` [friend]

## 27.166.5 Member Data Documentation

27.166.5.1 `int gdcM::JPEGCodec::BitSample` [protected]

27.166.5.2 `int gdcM::JPEGCodec::Quality` [protected]

The documentation for this class was generated from the following file:


- [gdcMJPEGCodec.h](#)

## 27.167 gdcM::JPEGLSCodec Class Reference


JPEG-LS.

```
#include <gdcmJPEGLSCodec.h>
```

Inheritance diagram for gdcm::JPEGLSCodec:

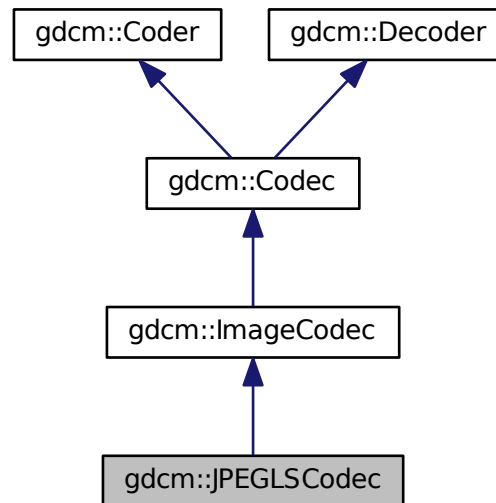

Collaboration diagram for gdcm::JPEGLSCodec:

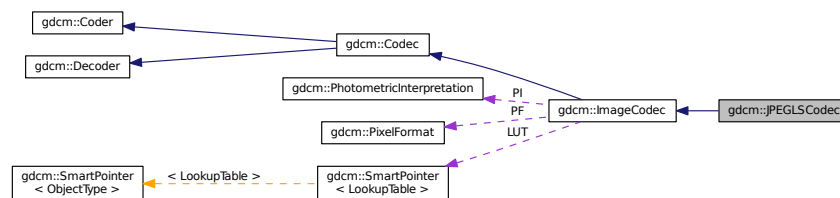

## Public Member Functions

- `JPEGLSCodec ()`
- `~JPEGLSCodec ()`
- `bool CanCode (TransferSyntax const &ts) const`  
*Return whether this coder support this transfer syntax (can code it)*
- `bool CanDecode (TransferSyntax const &ts) const`  
*Return whether this decoder support this transfer syntax (can decode it)*
- `virtual ImageCodec * Clone () const`
- `bool Code (DataElement const &in, DataElement &out)`  
*Code.*

- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)  
*Decode.*
  - bool [Decode](#) ([DataElement](#) const &in, char \*outBuffer, size\_t inBufferLength, uint32\_t inXMin, uint32\_t inXMax, uint32\_t inYMin, uint32\_t inYMax, uint32\_t inZMin, uint32\_t inZMax)
  - unsigned long [GetBufferLength](#) () const
  - bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
  - bool [GetLossless](#) () const
  - void [SetBufferLength](#) (unsigned long l)
  - void [SetLossless](#) (bool l)
  - void [SetLossyError](#) (int error)
- [0-3] generally*

### Protected Member Functions

- bool [AppendFrameEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [AppendRowEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [DecodeExtent](#) (char \*buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream &is)
- bool [IsFrameEncoder](#) ()
- bool [IsRowEncoder](#) ()
- bool [StartEncode](#) (std::ostream &)
- bool [StopEncode](#) (std::ostream &)

### Friends

- class [ImageRegionReader](#)

### Additional Inherited Members

#### 27.167.1 Detailed Description

JPEG-LS.

##### Note

codec that implement the JPEG-LS compression this is an implementation of [ImageCodec](#) for JPEG-LS

It uses the CharLS JPEG-LS implementation <http://charls.codeplex.com>

#### 27.167.2 Constructor & Destructor Documentation

27.167.2.1 `gdcm::JPEGLSCodec::JPEGLSCodec ( )`

27.167.2.2 `gdcm::JPEGLSCodec::~~JPEGLSCodec ( )`

#### 27.167.3 Member Function Documentation

27.167.3.1 `bool gdcm::JPEGLSCodec::AppendFrameEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.2 `bool gdcm::JPEGLSCodec::AppendRowEncode ( std::ostream & out, const char * data, size_t datalen )`  
`[protected], [virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.3 `bool gdcm::JPEGLSCodec::CanCode ( TransferSyntax const & ) const` `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.4 `bool gdcm::JPEGLSCodec::CanDecode ( TransferSyntax const & ) const` `[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.5 `virtual ImageCodec* gdcm::JPEGLSCodec::Clone ( ) const` `[virtual]`

Implements [gdcm::ImageCodec](#).

27.167.3.6 `bool gdcm::JPEGLSCodec::Code ( DataElement const & in_, DataElement & out_ )` `[virtual]`

Code.

Reimplemented from [gdcm::Coder](#).

27.167.3.7 `bool gdcm::JPEGLSCodec::Decode ( DataElement const & , DataElement & )` `[virtual]`

Decode.

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.8 `bool gdcm::JPEGLSCodec::Decode ( DataElement const & in, char * outBuffer, size_t inBufferLength, uint32_t inXMin, uint32_t inXMax, uint32_t inYMin, uint32_t inYMax, uint32_t inZMin, uint32_t inZMax )`

27.167.3.9 `bool gdcm::JPEGLSCodec::DecodeExtent ( char * buffer, unsigned int xmin, unsigned int xmax, unsigned int ymin, unsigned int ymax, unsigned int zmin, unsigned int zmax, std::istream & is )` `[protected]`

27.167.3.10 `unsigned long gdcm::JPEGLSCodec::GetBufferLength ( ) const` `[inline]`

27.167.3.11 `bool gdcm::JPEGLSCodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` `[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.12 `bool gdcm::JPEGLSCodec::GetLossless ( ) const`

27.167.3.13 `bool gdcm::JPEGLSCodec::IsFrameEncoder ( )` `[protected], [virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.167.3.14 `bool gdcml::JPEGLSCodec::IsRowEncoder ( )` [protected],[virtual]

Reimplemented from [gdcml::ImageCodec](#).

27.167.3.15 `void gdcml::JPEGLSCodec::SetBufferLength ( unsigned long / )` [inline]

27.167.3.16 `void gdcml::JPEGLSCodec::SetLossless ( bool / )`

27.167.3.17 `void gdcml::JPEGLSCodec::SetLossyError ( int error )`

[0-3] generally

27.167.3.18 `bool gdcml::JPEGLSCodec::StartEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcml::ImageCodec](#).

27.167.3.19 `bool gdcml::JPEGLSCodec::StopEncode ( std::ostream & )` [protected],[virtual]

Reimplemented from [gdcml::ImageCodec](#).

## 27.167.4 Friends And Related Function Documentation

27.167.4.1 `friend class ImageRegionReader` [friend]

The documentation for this class was generated from the following file:

- [gdcmlJPEGLSCodec.h](#)

## 27.168 gdcml::JSON Class Reference

```
#include <gdcmlJSON.h>
```

### Public Member Functions

- [JSON](#) ()
- [~JSON](#) ()
- bool [Code](#) ([DataSet](#) const &in, std::ostream &os)
- bool [Decode](#) (std::istream &is, [DataSet](#) &out)
- bool [GetPrettyPrint](#) () const
- void [PrettyPrintOff](#) ()
- void [PrettyPrintOn](#) ()
- void [SetPrettyPrint](#) (bool onoff)

### 27.168.1 Detailed Description

Examples:

[QIDO-RS.cxx](#).

## 27.168.2 Constructor & Destructor Documentation

27.168.2.1 `gdcm::JSON::JSON ( )`

27.168.2.2 `gdcm::JSON::~~JSON ( )`

## 27.168.3 Member Function Documentation

27.168.3.1 `bool gdcm::JSON::Code ( DataSet const & in, std::ostream & os )`

Examples:

[QIDO-RS.cxx](#).

27.168.3.2 `bool gdcm::JSON::Decode ( std::istream & is, DataSet & out )`

Examples:

[QIDO-RS.cxx](#).

27.168.3.3 `bool gdcm::JSON::GetPrettyPrint ( ) const`

27.168.3.4 `void gdcm::JSON::PrettyPrintOff ( )`

27.168.3.5 `void gdcm::JSON::PrettyPrintOn ( )`

Examples:

[QIDO-RS.cxx](#).

27.168.3.6 `void gdcm::JSON::SetPrettyPrint ( bool onoff )`

The documentation for this class was generated from the following file:


- [gdcmJSON.h](#)

## 27.169 gdcm::KAKADUCodec Class Reference


[KAKADUCodec](#).

```
#include <gdcmKAKADUCodec.h>
```

Inheritance diagram for `gdcm::KAKADUCodec`:

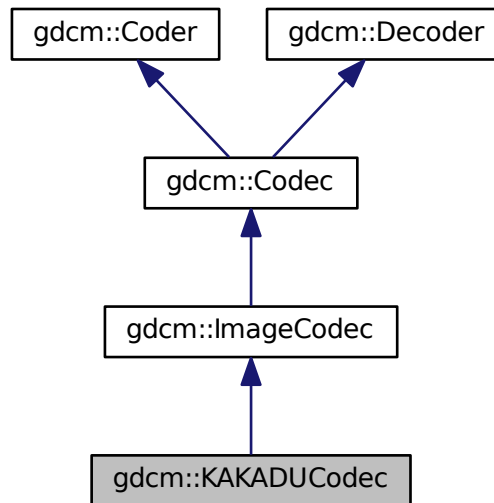

Collaboration diagram for `gdcm::KAKADUCodec`:

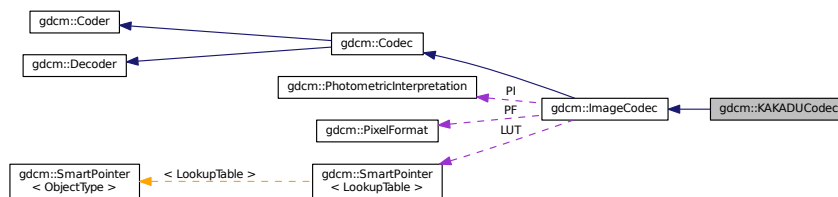

## Public Member Functions

- [`KAKADUCodec`](#) ()
- [`~KAKADUCodec`](#) ()
- [`bool CanCode`](#) ([`TransferSyntax`](#) const &ts) const  
Return whether this coder support this transfer syntax (can code it)
- [`bool CanDecode`](#) ([`TransferSyntax`](#) const &ts) const  
Return whether this decoder support this transfer syntax (can decode it)
- virtual [`ImageCodec \* Clone`](#) () const
- [`bool Code`](#) ([`DataElement`](#) const &in, [`DataElement`](#) &out)  
Code.
- [`bool Decode`](#) ([`DataElement`](#) const &is, [`DataElement`](#) &os)  
Decode.



## Additional Inherited Members

### 27.169.1 Detailed Description

[KAKADUCodec](#).

### 27.169.2 Constructor & Destructor Documentation

27.169.2.1 `gdcm::KAKADUCodec::KAKADUCodec ( )`

27.169.2.2 `gdcm::KAKADUCodec::~~KAKADUCodec ( )`

### 27.169.3 Member Function Documentation

27.169.3.1 `bool gdcm::KAKADUCodec::CanCode ( TransferSyntax const & ) const` `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.169.3.2 `bool gdcm::KAKADUCodec::CanDecode ( TransferSyntax const & ) const` `[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.169.3.3 `virtual ImageCodec* gdcm::KAKADUCodec::Clone ( ) const` `[virtual]`

Implements [gdcm::ImageCodec](#).

27.169.3.4 `bool gdcm::KAKADUCodec::Code ( DataElement const & in_, DataElement & out_ )` `[virtual]`

Code.

Reimplemented from [gdcm::Coder](#).

27.169.3.5 `bool gdcm::KAKADUCodec::Decode ( DataElement const &, DataElement & )` `[virtual]`

Decode.

Reimplemented from [gdcm::ImageCodec](#).

The documentation for this class was generated from the following file:


- [gdcmKAKADUCodec.h](#)

## 27.170 gdcm::LO Class Reference


[LO](#).

```
#include <gdcmLO.h>
```

Inheritance diagram for `gdcm::LO`:

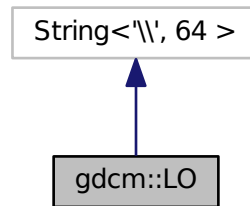

Collaboration diagram for `gdcm::LO`:

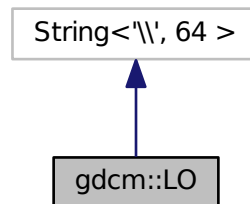

## Public Types

- typedef Superclass::const\_iterator [const\\_iterator](#)
- typedef Superclass::const\_reference [const\\_reference](#)
- typedef Superclass::const\_reverse\_iterator [const\\_reverse\\_iterator](#)
- typedef Superclass::difference\_type [difference\\_type](#)
- typedef Superclass::iterator [iterator](#)
- typedef Superclass::pointer [pointer](#)
- typedef Superclass::reference [reference](#)
- typedef Superclass::reverse\_iterator [reverse\\_iterator](#)
- typedef Superclass::size\_type [size\\_type](#)
- typedef [String<'\\', 64 >](#) [Superclass](#)
- typedef Superclass::value\_type [value\\_type](#)

## Public Member Functions

- [LO](#) ()

- [LO](#) (const [value\\_type](#) \*s)
- [LO](#) (const [value\\_type](#) \*s, [size\\_type](#) n)
- [LO](#) (const [Superclass](#) &s, [size\\_type](#) pos=0, [size\\_type](#) n=npos)
- bool [IsValid](#) () const

### 27.170.1 Detailed Description

[LO](#).

Note

TODO

### 27.170.2 Member Typedef Documentation

27.170.2.1 `typedef Superclass::const_iterator gdcm::LO::const_iterator`

27.170.2.2 `typedef Superclass::const_reference gdcm::LO::const_reference`

27.170.2.3 `typedef Superclass::const_reverse_iterator gdcm::LO::const_reverse_iterator`

27.170.2.4 `typedef Superclass::difference_type gdcm::LO::difference_type`

27.170.2.5 `typedef Superclass::iterator gdcm::LO::iterator`

27.170.2.6 `typedef Superclass::pointer gdcm::LO::pointer`

27.170.2.7 `typedef Superclass::reference gdcm::LO::reference`

27.170.2.8 `typedef Superclass::reverse_iterator gdcm::LO::reverse_iterator`

27.170.2.9 `typedef Superclass::size_type gdcm::LO::size_type`

27.170.2.10 `typedef String<'\',64> gdcm::LO::Superclass`

27.170.2.11 `typedef Superclass::value_type gdcm::LO::value_type`

### 27.170.3 Constructor & Destructor Documentation

27.170.3.1 `gdcm::LO::LO ( ) [inline]`

27.170.3.2 `gdcm::LO::LO ( const value_type * s ) [inline]`

27.170.3.3 `gdcm::LO::LO ( const value_type * s, size_type n ) [inline]`

27.170.3.4 `gdcm::LO::LO ( const Superclass & s, size_type pos = 0, size_type n = npos ) [inline]`

### 27.170.4 Member Function Documentation

27.170.4.1 `bool gdcm::LO::IsValid ( ) const [inline]`

The documentation for this class was generated from the following file:


- [gdcmLO.h](#)

## 27.171 gdcm::LookupTable Class Reference


[LookupTable](#) class.

```
#include <gdcmLookupTable.h>
```

Inheritance diagram for gdcm::LookupTable:

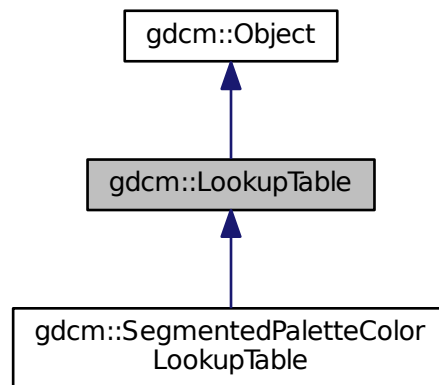

Collaboration diagram for gdcm::LookupTable:

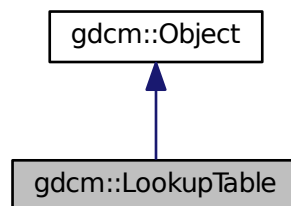

## Public Types

- enum [LookupTableType](#) {  
[RED](#) = 0,  
[GREEN](#),  
[BLUE](#),  
[GRAY](#),  
[UNKNOWN](#) }

## Public Member Functions

- [LookupTable](#) ()
- [LookupTable](#) ([LookupTable](#) const &lut)
- [~LookupTable](#) ()
- void [Allocate](#) (unsigned short bitsample=8)  
*Allocate the LUT.*
- void [Clear](#) ()  
*Clear the LUT.*
- void [Decode](#) (std::istream &is, std::ostream &os) const  
*Decode the LUT.*
- bool [Decode](#) (char \*outputbuffer, size\_t outlen, const char \*inputbuffer, size\_t inlen) const
- unsigned short [GetBitSample](#) () const  
*return the bit sample*
- bool [GetBufferAsRGBA](#) (unsigned char \*rgba) const  
*return the LUT as RGBA buffer*
- void [GetLUT](#) ([LookupTableType](#) type, unsigned char \*array, unsigned int &length) const
- void [GetLUTDescriptor](#) ([LookupTableType](#) type, unsigned short &length, unsigned short &subscript, unsigned short &bitsize) const
- unsigned int [GetLUTLength](#) ([LookupTableType](#) type) const
- const unsigned char \* [GetPointer](#) () const  
*return a raw pointer to the LUT*
- void [InitializeBlueLUT](#) (unsigned short length, unsigned short subscript, unsigned short bitsize)
- bool [Initialized](#) () const  
*return whether the LUT has been initialized*
- void [InitializeGreenLUT](#) (unsigned short length, unsigned short subscript, unsigned short bitsize)
- void [InitializeLUT](#) ([LookupTableType](#) type, unsigned short length, unsigned short subscript, unsigned short bitsize)  
*Generic interface:*
- void [InitializeRedLUT](#) (unsigned short length, unsigned short subscript, unsigned short bitsize)  
*RED / GREEN / BLUE specific:*
- void [Print](#) (std::ostream &) const
- void [SetBlueLUT](#) (const unsigned char \*blue, unsigned int length)
- void [SetGreenLUT](#) (const unsigned char \*green, unsigned int length)
- virtual void [SetLUT](#) ([LookupTableType](#) type, const unsigned char \*array, unsigned int length)
- void [SetRedLUT](#) (const unsigned char \*red, unsigned int length)
- bool [WriteBufferAsRGBA](#) (const unsigned char \*rgba)  
*Write the LUT as RGBA.*

## Protected Attributes

- unsigned short [BitSample](#)
- bool [IncompleteLUT](#):1
- LookupTableInternal \* [Internal](#)

## Additional Inherited Members

### 27.171.1 Detailed Description

[LookupTable](#) class.

### 27.171.2 Member Enumeration Documentation

#### 27.171.2.1 enum gdcm::LookupTable::LookupTableType

Enumerator

***RED***

***GREEN***

***BLUE***

***GRAY***

***UNKNOWN***

### 27.171.3 Constructor & Destructor Documentation

27.171.3.1 gdcm::LookupTable::LookupTable ( )

27.171.3.2 gdcm::LookupTable::~~LookupTable ( )

27.171.3.3 gdcm::LookupTable::LookupTable ( [LookupTable](#) const & *lut* ) `[inline]`

### 27.171.4 Member Function Documentation

27.171.4.1 void gdcm::LookupTable::Allocate ( unsigned short *bitsample* = 8 )

Allocate the LUT.

27.171.4.2 void gdcm::LookupTable::Clear ( )

Clear the LUT.

27.171.4.3 void gdcm::LookupTable::Decode ( std::istream & *is*, std::ostream & *os* ) const

Decode the LUT.

27.171.4.4 `bool gdcm::LookupTable::Decode ( char * outputbuffer, size_t outlen, const char * inputbuffer, size_t inlen ) const`

Decode the LUT outputbuffer will contains the RGB decoded PALETTE COLOR input image of size inlen the outputbuffer should be at least 3 times the size of inlen

27.171.4.5 `unsigned short gdcm::LookupTable::GetBitSample ( ) const [inline]`

return the bit sample

27.171.4.6 `bool gdcm::LookupTable::GetBufferAsRGBA ( unsigned char * rgba ) const`

return the LUT as RGBA buffer

27.171.4.7 `void gdcm::LookupTable::GetLUT ( LookupTableType type, unsigned char * array, unsigned int & length ) const`

Examples:

[ExtractImageRegionWithLUT.cs](#).

27.171.4.8 `void gdcm::LookupTable::GetLUTDescriptor ( LookupTableType type, unsigned short & length, unsigned short & subscript, unsigned short & bitsize ) const`

27.171.4.9 `unsigned int gdcm::LookupTable::GetLUTLength ( LookupTableType type ) const`

27.171.4.10 `const unsigned char* gdcm::LookupTable::GetPointer ( ) const`

return a raw pointer to the LUT

27.171.4.11 `void gdcm::LookupTable::InitializeBlueLUT ( unsigned short length, unsigned short subscript, unsigned short bitsize )`

27.171.4.12 `bool gdcm::LookupTable::Initialized ( ) const`

return whether the LUT has been initialized

27.171.4.13 `void gdcm::LookupTable::InitializeGreenLUT ( unsigned short length, unsigned short subscript, unsigned short bitsize )`

27.171.4.14 `void gdcm::LookupTable::InitializeLUT ( LookupTableType type, unsigned short length, unsigned short subscript, unsigned short bitsize )`

Generic interface:

27.171.4.15 `void gdcm::LookupTable::InitializeRedLUT ( unsigned short length, unsigned short subscript, unsigned short bitsize )`

RED / GREEN / BLUE specific:

27.171.4.16 `void gdcm::LookupTable::Print ( std::ostream & ) const` `[inline],[virtual]`

Reimplemented from [gdcm::Object](#).

Reimplemented in [gdcm::SegmentedPaletteColorLookupTable](#).

27.171.4.17 `void gdcm::LookupTable::SetBlueLUT ( const unsigned char * blue, unsigned int length )`

27.171.4.18 `void gdcm::LookupTable::SetGreenLUT ( const unsigned char * green, unsigned int length )`

27.171.4.19 `virtual void gdcm::LookupTable::SetLUT ( LookupTableType type, const unsigned char * array, unsigned int length )` `[virtual]`

Reimplemented in [gdcm::SegmentedPaletteColorLookupTable](#).

27.171.4.20 `void gdcm::LookupTable::SetRedLUT ( const unsigned char * red, unsigned int length )`

27.171.4.21 `bool gdcm::LookupTable::WriteBufferAsRGBA ( const unsigned char * rgba )`

Write the LUT as RGBA.

## 27.171.5 Member Data Documentation

27.171.5.1 `unsigned short gdcm::LookupTable::BitSample` `[protected]`

27.171.5.2 `bool gdcm::LookupTable::IncompleteLUT` `[protected]`

27.171.5.3 `LookupTableInternal* gdcm::LookupTable::Internal` `[protected]`

The documentation for this class was generated from the following file:

- [gdcmLookupTable.h](#)

## 27.172 gdcm::Scanner::Itstr Struct Reference

```
#include <gdcmScanner.h>
```

### Public Member Functions

- `bool operator() (const char *s1, const char *s2) const`

## 27.172.1 Member Function Documentation

27.172.1.1 `bool gdcm::Scanner::Itstr::operator() ( const char * s1, const char * s2 ) const` `[inline]`

The documentation for this struct was generated from the following file:

- [gdcmScanner.h](#)



## 27.173 gdcmmacro::Macro Class Reference

Class for representing a [Macro](#).

```
#include <gdcmmacro.h>
```

### Public Types

- typedef std::vector< std::string > [ArrayIncludeMacrosType](#)
- typedef std::map< [Tag](#), [MacroEntry](#) > [MapModuleEntry](#)

### Public Member Functions

- [Macro](#) ()
- void [AddMacroEntry](#) (const [Tag](#) &tag, const [MacroEntry](#) &module)  
*Will add a [ModuleEntry](#) directly at root-level. See [Macro](#) for nested-included level.*
- void [Clear](#) ()
- bool [FindMacroEntry](#) (const [Tag](#) &tag) const
- const [MacroEntry](#) & [GetMacroEntry](#) (const [Tag](#) &tag) const
- const char \* [GetName](#) () const
- void [SetName](#) (const char \*name)
- bool [Verify](#) (const [DataSet](#) &ds, [Usage](#) const &usage) const

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Macro](#) &\_val)

#### 27.173.1 Detailed Description

Class for representing a [Macro](#).

#### Note

**Attribute Macro:** a set of Attributes that are described in a single table that is referenced by multiple [Module](#) or other tables.

#### See also

[Module](#)

#### 27.173.2 Member Typedef Documentation

27.173.2.1 typedef std::vector<std::string> gdcmmacro::Macro::ArrayIncludeMacrosType

27.173.2.2 typedef std::map<Tag, MacroEntry> gdcmmacro::Macro::MapModuleEntry

#### 27.173.3 Constructor & Destructor Documentation

27.173.3.1 gdcmmacro::Macro ( ) [inline]

### 27.173.4 Member Function Documentation

27.173.4.1 `void gdcmmacro::Macro::AddMacroEntry ( const Tag & tag, const MacroEntry & module ) [inline]`

Will add a [ModuleEntry](#) directly at root-level. See [Macro](#) for nested-included level.

27.173.4.2 `void gdcmmacro::Macro::Clear ( ) [inline]`

27.173.4.3 `bool gdcmmacro::Macro::FindMacroEntry ( const Tag & tag ) const`

Find or Get a [ModuleEntry](#). [ModuleEntry](#) are either search are root-level or within nested-macro included in module.

27.173.4.4 `const MacroEntry& gdcmmacro::Macro::GetMacroEntry ( const Tag & tag ) const`

27.173.4.5 `const char* gdcmmacro::Macro::GetName ( ) const [inline]`

27.173.4.6 `void gdcmmacro::Macro::SetName ( const char * name ) [inline]`

27.173.4.7 `bool gdcmmacro::Macro::Verify ( const DataSet & ds, Usage const & usage ) const`

### 27.173.5 Friends And Related Function Documentation

27.173.5.1 `std::ostream& operator<< ( std::ostream & _os, const Macro & _val ) [friend]`

The documentation for this class was generated from the following file:

- [gdcmmacro.h](#)

## 27.174 gdcmmacro::Macros Class Reference

Class for representing a [Modules](#).

```
#include <gdcmmacros.h>
```

### Public Types

- `typedef std::map< std::string, Macro > ModuleMapType`

### Public Member Functions

- [Macros](#) ()
- `void AddMacro (const char *ref, const Macro &module)`
- `void Clear ()`
- `const Macro & GetMacro (const char *name) const`
- `bool IsEmpty () const`

### Friends

- `std::ostream & operator<< (std::ostream & _os, const Macros & _val)`

### 27.174.1 Detailed Description

Class for representing a [Modules](#).

#### Note

bla

#### See also

[Module](#)

#### Examples:

[TraverseModules.cxx](#).

### 27.174.2 Member Typedef Documentation

27.174.2.1 `typedef std::map<std::string, Macro> gdcm::Macros::ModuleMapType`

### 27.174.3 Constructor & Destructor Documentation

27.174.3.1 `gdcm::Macros::Macros ( )` `[inline]`

### 27.174.4 Member Function Documentation

27.174.4.1 `void gdcm::Macros::AddMacro ( const char * ref, const Macro & module )` `[inline]`

27.174.4.2 `void gdcm::Macros::Clear ( )` `[inline]`

27.174.4.3 `const Macro& gdcm::Macros::GetMacro ( const char * name ) const` `[inline]`

27.174.4.4 `bool gdcm::Macros::IsEmpty ( ) const` `[inline]`

### 27.174.5 Friends And Related Function Documentation

27.174.5.1 `std::ostream& operator<< ( std::ostream & _os, const Macros & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmMacros.h](#)

## 27.175 gdcm::network::MaximumLengthSub Class Reference

[MaximumLengthSub](#) Annex D [Table D.1-1](#) MAXIMUM LENGTH SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

```
#include <gdcmMaximumLengthSub.h>
```

### Public Member Functions

- [MaximumLengthSub](#) ()

- uint32\_t [GetMaximumLength](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetMaximumLength](#) (uint32\_t maximumlength)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.175.1 Detailed Description

[MaximumLengthSub](#) Annex D [Table D.1-1](#) MAXIMUM LENGTH SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

or

[Table D.1-2](#) Maximum length sub-item fields (A-ASSOCIATE-AC)

### 27.175.2 Constructor & Destructor Documentation

27.175.2.1 `gdcmm::network::MaximumLengthSub::MaximumLengthSub ( )`

### 27.175.3 Member Function Documentation

27.175.3.1 `uint32_t gdcmm::network::MaximumLengthSub::GetMaximumLength ( ) const` `[inline]`

27.175.3.2 `void gdcmm::network::MaximumLengthSub::Print ( std::ostream & os ) const`

27.175.3.3 `std::istream& gdcmm::network::MaximumLengthSub::Read ( std::istream & is )`

27.175.3.4 `void gdcmm::network::MaximumLengthSub::SetMaximumLength ( uint32_t maximumlength )`

27.175.3.5 `size_t gdcmm::network::MaximumLengthSub::Size ( ) const`

27.175.3.6 `const std::ostream& gdcmm::network::MaximumLengthSub::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:

- [gdcmmMaximumLengthSub.h](#)

## 27.176 gdcmm::MD5 Class Reference

Class for [MD5](#).

```
#include <gdcmmMD5.h>
```

### Public Member Functions

- [MD5](#) ()
- [~MD5](#) ()

## Static Public Member Functions

- static bool [Compute](#) (const char \*buffer, unsigned long buf\_len, char digest\_str[33])
- static bool [ComputeFile](#) (const char \*filename, char digest\_str[33])

### 27.176.1 Detailed Description

Class for [MD5](#).

#### Warning

this class is able to pick from two implementations:

1. a lightweight md5 implementation (when GDCM\_BUILD\_TESTING is turned ON)
2. the one from OpenSSL (when GDCM\_USE\_SYSTEM\_OPENSSL is turned ON)

In all other cases it will return an error

### 27.176.2 Constructor & Destructor Documentation

27.176.2.1 `gdcm::MD5::MD5 ( )`

27.176.2.2 `gdcm::MD5::~~MD5 ( )`

### 27.176.3 Member Function Documentation

27.176.3.1 `static bool gdcm::MD5::Compute ( const char * buffer, unsigned long buf_len, char digest_str[33] )` `[static]`

27.176.3.2 `static bool gdcm::MD5::ComputeFile ( const char * filename, char digest_str[33] )` `[static]`

The documentation for this class was generated from the following file:

- [gdcmMD5.h](#)

## 27.177 gdcm::MediaStorage Class Reference

[MediaStorage](#).

```
#include <gdcmMediaStorage.h>
```

## Public Types

- enum `MSType` {
  - `MediaStorageDirectoryStorage` = 0,
  - `ComputedRadiographyImageStorage`,
  - `DigitalXRayImageStorageForPresentation`,
  - `DigitalXRayImageStorageForProcessing`,
  - `DigitalMammographyImageStorageForPresentation`,
  - `DigitalMammographyImageStorageForProcessing`,
  - `DigitalIntraoralXRayImageStorageForPresentation`,
  - `DigitalIntraoralXRayImageStorageForProcessing`,
  - `CTImageStorage`,
  - `EnhancedCTImageStorage`,
  - `UltrasoundImageStorageRetired`,
  - `UltrasoundImageStorage`,
  - `UltrasoundMultiFrameImageStorageRetired`,
  - `UltrasoundMultiFrameImageStorage`,
  - `MRImageStorage`,
  - `EnhancedMRImageStorage`,
  - `MRSpectroscopyStorage`,
  - `NuclearMedicineImageStorageRetired`,
  - `SecondaryCaptureImageStorage`,
  - `MultiframeSingleBitSecondaryCaptureImageStorage`,
  - `MultiframeGrayscaleByteSecondaryCaptureImageStorage`,
  - `MultiframeGrayscaleWordSecondaryCaptureImageStorage`,
  - `MultiframeTrueColorSecondaryCaptureImageStorage`,
  - `StandaloneOverlayStorage`,
  - `StandaloneCurveStorage`,
  - `LeadECGWaveformStorage`,
  - `GeneralECGWaveformStorage`,
  - `AmbulatoryECGWaveformStorage`,
  - `HemodynamicWaveformStorage`,
  - `CardiacElectrophysiologyWaveformStorage`,
  - `BasicVoiceAudioWaveformStorage`,
  - `StandaloneModalityLUTStorage`,
  - `StandaloneVOILUTStorage`,
  - `GrayscaleSoftcopyPresentationStateStorageSOPClass`,
  - `XRayAngiographicImageStorage`,
  - `XRayRadiofluoroscopicImageStorage`,
  - `XRayAngiographicBiPlaneImageStorageRetired`,
  - `NuclearMedicineImageStorage`,
  - `RawDataStorage`,
  - `SpacialRegistrationStorage`,
  - `SpacialFiducialsStorage`,
  - `PETImageStorage`,
  - `RTImageStorage`,
  - `RTDoseStorage`,
  - `RTStructureSetStorage`,
  - `RTPlanStorage`,
  - `CSANonImageStorage`,
  - `Philips3D`,
  - `EnhancedSR`,
  - `BasicTextSR`,
  - `HardcopyGrayscaleImageStorage`,
  - `ComprehensiveSR`,
  - `DetachedStudyManagementSOPClass`,
  - `EncapsulatedPDFStorage`,
  - `EncapsulatedCDASStorage`,
  - `StudyComponentManagementSOPClass`,
  - `DetachedVisitManagementSOPClass`,
  - `DetachedPatientManagementSOPClass`,

MS\_END }

- enum `ObjectType` {  
`NoObject` = 0,  
`Video`,  
`Waveform`,  
`Audio`,  
`PDF`,  
`URI`,  
`Segmentation`,  
`ObjectEnd` }

## Public Member Functions

- `MediaStorage` (`MSType` type=`MS_END`)
- const char \* `GetModality` () const
- unsigned int `GetModalityDimension` () const
- const char \* `GetString` () const  
*Return the Media [String](#) of the object.*
- void `GuessFromModality` (const char \*modality, unsigned int dimension=2)
- bool `IsUndefined` () const
- operator `MSType` () const
- bool `SetFromDataSet` (`DataSet` const &ds)
- bool `SetFromFile` (`File` const &file)
- bool `SetFromHeader` (`FileMetaInformation` const &fmi)
- bool `SetFromModality` (`DataSet` const &ds)

## Static Public Member Functions

- static const char \* `GetMSString` (`MSType` ts)  
*Return the Media [String](#) associated. Will return NULL for MS\_END.*
- static `MSType` `GetMSType` (const char \*str)
- static unsigned int `GetNumberOfModality` ()
- static unsigned int `GetNumberOfMSString` ()
- static unsigned int `GetNumberOfMSType` ()
- static bool `IsImage` (`MSType` ts)

## Protected Member Functions

- void `SetFromSourceImageSequence` (`DataSet` const &ds)

## Friends

- std::ostream & `operator<<` (std::ostream &os, const `MediaStorage` &ms)

### 27.177.1 Detailed Description

[MediaStorage](#).

#### Note

FIXME There should not be any notion of [Image](#) and/or PDF at that point Only the codec can answer yes I support this Media Storage or not... For instance an [ImageCodec](#) will answer yes to most of them while a [PDFCodec](#) will answer only for the Encapsulated PDF

#### See also

[UIDs](#)

#### Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllVR.cxx](#), [GenerateStandardSOPClasses.cxx](#), [GenFakeIdentifyFile.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [StreamImageReaderTest.cxx](#), and [TestReader.cxx](#).

### 27.177.2 Member Enumeration Documentation

#### 27.177.2.1 enum gdcm::MediaStorage::MSType

##### Enumerator

***MediaStorageDirectoryStorage***  
***ComputedRadiographyImageStorage***  
***DigitalXRayImageStorageForPresentation***  
***DigitalXRayImageStorageForProcessing***  
***DigitalMammographyImageStorageForPresentation***  
***DigitalMammographyImageStorageForProcessing***  
***DigitalIntraoralXrayImageStorageForPresentation***  
***DigitalIntraoralXRayImageStorageForProcessing***  
***CTImageStorage***  
***EnhancedCTImageStorage***  
***UltrasoundImageStorageRetired***  
***UltrasoundImageStorage***  
***UltrasoundMultiFrameImageStorageRetired***  
***UltrasoundMultiFrameImageStorage***  
***MRImageStorage***  
***EnhancedMRImageStorage***  
***MRSpectroscopyStorage***  
***NuclearMedicineImageStorageRetired***  
***SecondaryCaptureImageStorage***  
***MultiframeSingleBitSecondaryCaptureImageStorage***  
***MultiframeGrayscaleByteSecondaryCaptureImageStorage***



*MultiframeGrayscaleWordSecondaryCaptureImageStorage*  
*MultiframeTrueColorSecondaryCaptureImageStorage*  
*StandaloneOverlayStorage*  
*StandaloneCurveStorage*  
*LeadECGWaveformStorage*  
*GeneralECGWaveformStorage*  
*AmbulatoryECGWaveformStorage*  
*HemodynamicWaveformStorage*  
*CardiacElectrophysiologyWaveformStorage*  
*BasicVoiceAudioWaveformStorage*  
*StandaloneModalityLUTStorage*  
*StandaloneVOILUTStorage*  
*GrayscaleSoftcopyPresentationStateStorageSOPClass*  
*XRayAngiographicImageStorage*  
*XRayRadiofluoroscopicImageStorage*  
*XRayAngiographicBiPlaneImageStorageRetired*  
*NuclearMedicineImageStorage*  
*RawDataStorage*  
*SpacialRegistrationStorage*  
*SpacialFiducialsStorage*  
*PETImageStorage*  
*RTImageStorage*  
*RTDoseStorage*  
*RTStructureSetStorage*  
*RTPlanStorage*  
*CSANonImageStorage*  
*Philips3D*  
*EnhancedSR*  
*BasicTextSR*  
*HardcopyGrayscaleImageStorage*  
*ComprehensiveSR*  
*DetachedStudyManagementSOPClass*  
*EncapsulatedPDFStorage*  
*EncapsulatedCDASStorage*  
*StudyComponentManagementSOPClass*  
*DetachedVisitManagementSOPClass*  
*DetachedPatientManagementSOPClass*  
*VideoEndoscopicImageStorage*  
*GeneralElectricMagneticResonanceImageStorage*  
*GEPrivate3DModelStorage*  
*ToshibaPrivateDataStorage*  
*MammographyCADSR*

*KeyObjectSelectionDocument*  
*HangingProtocolStorage*  
*ModalityPerformedProcedureStepSOPClass*  
*PhilipsPrivateMRSyntheticImageStorage*  
*VLPhotographicImageStorage*  
*SegmentationStorage*  
*RTIonPlanStorage*  
*XR3DAngiographicImageStorage*  
*EnhancedXAImageStorage*  
*RTIonBeamsTreatmentRecordStorage*  
*SurfaceSegmentationStorage*  
*VLWholeSlideMicroscopyImageStorage*  
*RTTreatmentSummaryRecordStorage*  
*EnhancedUSVolumeStorage*  
*XR3DRadiationDoseSR*  
*VLEndoscopicImageStorage*  
*BreastTomosynthesisImageStorage*  
*FujiPrivateCRLImageStorage*  
*OphthalmicPhotography8BitImageStorage*  
*OphthalmicTomographyImageStorage*  
*VLMicroscopicImageStorage*  
*MS\_END*

Examples:

[GenerateStandardSOPClasses.cxx](#).

27.177.2.2 enum gdcmm::MediaStorage::ObjectType

Enumerator

*NoObject*  
*Video*  
*Waveform*  
*Audio*  
*PDF*  
*URI*  
*Segmentation*  
*ObjectEnd*

### 27.177.3 Constructor & Destructor Documentation

27.177.3.1 `gdcm::MediaStorage::MediaStorage ( MStype type = MS_END ) [inline]`

### 27.177.4 Member Function Documentation

27.177.4.1 `const char* gdcm::MediaStorage::GetModality ( ) const`

27.177.4.2 `unsigned int gdcm::MediaStorage::GetModalityDimension ( ) const`

27.177.4.3 `static const char* gdcm::MediaStorage::GetMSString ( MStype ts ) [static]`

Return the Media [String](#) associated. Will return NULL for MS\_END.

Examples:

[GenerateStandardSOPClasses.cxx](#).

Referenced by `gdcm::operator<<()`.

27.177.4.4 `static MStype gdcm::MediaStorage::GetMStype ( const char * str ) [static]`

Examples:

[MetaImageMD5Activiz.cs](#), and [TestReader.cxx](#).

27.177.4.5 `static unsigned int gdcm::MediaStorage::GetNumberOfModality ( ) [static]`

27.177.4.6 `static unsigned int gdcm::MediaStorage::GetNumberOfMSString ( ) [static]`

27.177.4.7 `static unsigned int gdcm::MediaStorage::GetNumberOfMStype ( ) [static]`

27.177.4.8 `const char* gdcm::MediaStorage::GetString ( ) const`

Return the Media [String](#) of the object.

Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), and [StreamImageReaderTest.cxx](#).

27.177.4.9 `void gdcm::MediaStorage::GuessFromModality ( const char * modality, unsigned int dimension = 2 )`

27.177.4.10 `static bool gdcm::MediaStorage::IsImage ( MStype ts ) [static]`

Returns whether DICOM has a Pixel Data element (7fe0,0010)

**Warning**

MRSpectroscopyStorage could be image but are not

**Examples:**

[MetalImageMD5Activiz.cs](#).

27.177.4.11 `bool gdcM::MediaStorage::IsUndefined ( ) const [inline]`

**Examples:**

[TestReader.cxx](#).

27.177.4.12 `gdcM::MediaStorage::operator MStype ( ) const [inline]`

27.177.4.13 `bool gdcM::MediaStorage::SetFromDataSet ( DataSet const & ds )`

Advanced user only (functions should be protected level...) Those function are lower level than SetFromFile

27.177.4.14 `bool gdcM::MediaStorage::SetFromFile ( File const & file )`

Attempt to set the [MediaStorage](#) from a file: WARNING: When no [MediaStorage](#) & Modality are found BUT a PixelData element is found then [MediaStorage](#) is set to the default SecondaryCaptureImageStorage (return value is false in this case)

**Examples:**

[gdcMrtionplan.cxx](#), [gdcMrtplan.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), and [TestReader.cxx](#).

27.177.4.15 `bool gdcM::MediaStorage::SetFromHeader ( FileMetaInformation const & fmi )`

27.177.4.16 `bool gdcM::MediaStorage::SetFromModality ( DataSet const & ds )`

27.177.4.17 `void gdcM::MediaStorage::SetFromSourceImageSequence ( DataSet const & ds ) [protected]`

**27.177.5 Friends And Related Function Documentation**

27.177.5.1 `std::ostream& operator<< ( std::ostream & os, const MediaStorage & ms ) [friend]`

The documentation for this class was generated from the following file:


- [gdcMMediaStorage.h](#)

**27.178 gdcM::MemberCommand< T > Class Template Reference**


[Command](#) subclass that calls a pointer to a member function.

```
#include <gdcMCommand.h>
```

Inheritance diagram for gdcM::MemberCommand< T >:

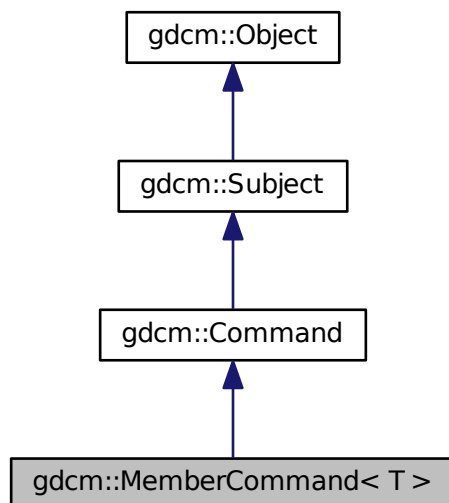

Collaboration diagram for gdcM::MemberCommand< T >:

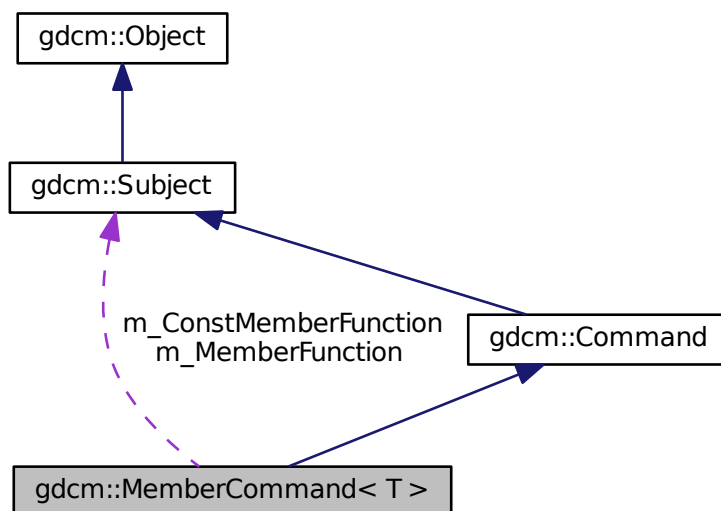

## Public Types

- typedef [MemberCommand](#) Self
- typedef void(T::\* [TConstMemberFunctionPointer](#)) (const [Subject](#) \*, const [Event](#) &)
- typedef void(T::\* [TMemberFunctionPointer](#)) ([Subject](#) \*, const [Event](#) &)

## Public Member Functions

- virtual void [Execute](#) ([Subject](#) \*caller, const [Event](#) &event)
- virtual void [Execute](#) (const [Subject](#) \*caller, const [Event](#) &event)
- void [SetCallbackFunction](#) (T \*object, [TMemberFunctionPointer](#) memberFunction)
- void [SetCallbackFunction](#) (T \*object, [TConstMemberFunctionPointer](#) memberFunction)

## Static Public Member Functions

- static [SmartPointer](#)< [MemberCommand](#) > [New](#) ()

## Protected Member Functions

- [MemberCommand](#) ()
- virtual [~MemberCommand](#) ()

## Protected Attributes

- [TConstMemberFunctionPointer](#) m\_ConstMemberFunction
- [TMemberFunctionPointer](#) m\_MemberFunction
- T \* m\_This

### 27.178.1 Detailed Description

template<class T>class [gdcmmembercommand](#)< T >

[Command](#) subclass that calls a pointer to a member function.

[MemberCommand](#) calls a pointer to a member function with the same arguments as [Execute](#) on [Command](#).

### 27.178.2 Member Typedef Documentation

27.178.2.1 template<class T > typedef [MemberCommand](#) [gdcmmembercommand](#)< T >::Self

Standard class typedefs.

27.178.2.2 template<class T > typedef void(T::\* [gdcmmembercommand](#)< T >::TConstMemberFunctionPointer) (const [Subject](#) \*, const [Event](#) &)

27.178.2.3 template<class T > typedef void(T::\* [gdcmmembercommand](#)< T >::TMemberFunctionPointer) ([Subject](#) \*, const [Event](#) &)

pointer to a member function that takes a [Subject](#)\* and the event

### 27.178.3 Constructor & Destructor Documentation

27.178.3.1 `template<class T> gdcmmembercommand< T >::membercommand ( ) [inline],  
[protected]`

Referenced by `gdcmmembercommand< T >::New()`.

27.178.3.2 `template<class T> virtual gdcmmembercommand< T >::~~membercommand ( ) [inline],  
[protected], [virtual]`

### 27.178.4 Member Function Documentation

27.178.4.1 `template<class T> virtual void gdcmmembercommand< T >::Execute ( Subject * caller, const Event &  
event ) [inline], [virtual]`

Invoke the member function.

Implements [gdcmmembercommand::Command](#).

References `gdcmmembercommand< T >::m_memberFunction`.

27.178.4.2 `template<class T> virtual void gdcmmembercommand< T >::Execute ( const Subject * caller, const Event  
& event ) [inline], [virtual]`

Invoke the member function with a const object.

Implements [gdcmmembercommand::Command](#).

References `gdcmmembercommand< T >::m_ConstMemberFunction`.

27.178.4.3 `template<class T> static SmartPointer<membercommand> gdcmmembercommand< T >::New ( )  
[inline], [static]`

Method for creation through the object factory.

References `gdcmmembercommand< T >::MemberCommand()`.

27.178.4.4 `template<class T> void gdcmmembercommand< T >::SetCallbackFunction ( T * object,  
TMemberFunctionPointer memberFunction ) [inline]`

Run-time type information (and related methods). Set the callback function along with the object that it will be invoked on.

References `gdcmmembercommand< T >::m_memberFunction`, and `gdcmmembercommand< T >::m_This`.

27.178.4.5 `template<class T> void gdcmmembercommand< T >::SetCallbackFunction ( T * object,  
TConstMemberFunctionPointer memberFunction ) [inline]`

References `gdcmmembercommand< T >::m_ConstMemberFunction`, and `gdcmmembercommand< T >::m_This`.

### 27.178.5 Member Data Documentation

27.178.5.1 `template<class T> TConstMemberFunctionPointer gdcM::MemberCommand< T  
>::m_ConstMemberFunction [protected]`

Referenced by `gdcM::MemberCommand< T >::Execute()`, and `gdcM::MemberCommand< T >::SetCallback←  
Function()`.

27.178.5.2 `template<class T> TMemberFunctionPointer gdcM::MemberCommand< T >::m_MemberFunction  
[protected]`

Referenced by `gdcM::MemberCommand< T >::Execute()`, and `gdcM::MemberCommand< T >::SetCallback←  
Function()`.

27.178.5.3 `template<class T> T* gdcM::MemberCommand< T >::m_This [protected]`

Referenced by `gdcM::MemberCommand< T >::SetCallbackFunction()`.

The documentation for this class was generated from the following file:


- [gdcMCommand.h](#)

## 27.179 gdcM::MeshPrimitive Class Reference

This class defines surface mesh primitives. It is designed from surface mesh primitives macro.


```
#include <gdcMMeshPrimitive.h>
```

Inheritance diagram for `gdcM::MeshPrimitive`:

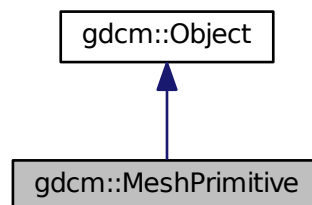



Collaboration diagram for gdcm::MeshPrimitive:

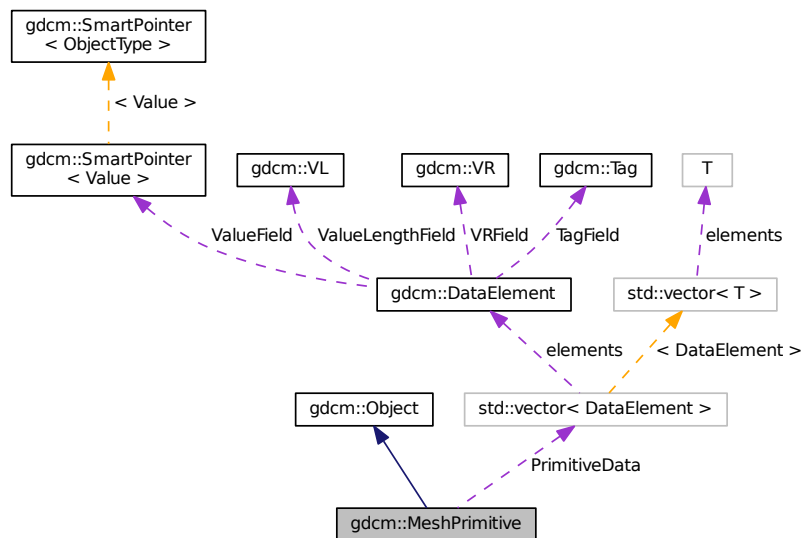

## Public Types

- enum `MPTType` {  
`VERTEX` = 0,  
`EDGE`,  
`TRIANGLE`,  
`TRIANGLE_STRIP`,  
`TRIANGLE_FAN`,  
`LINE`,  
`FACET`,  
`MPTType_END` }
- *This enumeration defines primitive types.*
- typedef `std::vector< DataElement >` `PrimitivesData`

## Public Member Functions

- `MeshPrimitive` ()
- virtual `~MeshPrimitive` ()
- void `AddPrimitiveData` (`DataElement` const &de)
- unsigned int `GetNumberOfPrimitivesData` () const
- const `DataElement` & `GetPrimitiveData` () const
- `DataElement` & `GetPrimitiveData` ()
- const `DataElement` & `GetPrimitiveData` (const unsigned int idx) const
- `DataElement` & `GetPrimitiveData` (const unsigned int idx)
- const `PrimitivesData` & `GetPrimitivesData` () const
- `PrimitivesData` & `GetPrimitivesData` ()
- `MPTType` `GetPrimitiveType` () const

- void [SetPrimitiveData](#) ([DataElement](#) const &de)
- void [SetPrimitiveData](#) (const unsigned int idx, [DataElement](#) const &de)
- void [SetPrimitivesData](#) ([PrimitivesData](#) const &DEs)
- void [SetPrimitiveType](#) (const [MPTType](#) type)

### Static Public Member Functions

- static [MPTType](#) [GetMPTType](#) (const char \*type)
- static const char \* [GetMPTTypeString](#) (const [MPTType](#) type)

### Protected Attributes

- [PrimitivesData](#) [PrimitiveData](#)
- [MPTType](#) [PrimitiveType](#)

### Additional Inherited Members

#### 27.179.1 Detailed Description

This class defines surface mesh primitives. It is designed from surface mesh primitives macro.

See also

PS 3.3 C.27.4

#### 27.179.2 Member Typedef Documentation

27.179.2.1 `typedef std::vector< DataElement > gdcm::MeshPrimitive::PrimitivesData`

#### 27.179.3 Member Enumeration Documentation

27.179.3.1 `enum gdcm::MeshPrimitive::MPTType`

This enumeration defines primitive types.

See also

PS 3.3 C.27.4.1

Enumerator

***VERTEX***

***EDGE***

***TRIANGLE***

***TRIANGLE\_STRIP***

***TRIANGLE\_FAN***

***LINE***

***FACET***

***MPTType\_END***

### 27.179.4 Constructor & Destructor Documentation

27.179.4.1 `gdcm::MeshPrimitive::MeshPrimitive ( )`

27.179.4.2 `virtual gdcm::MeshPrimitive::~~MeshPrimitive ( )` `[virtual]`

### 27.179.5 Member Function Documentation

27.179.5.1 `void gdcm::MeshPrimitive::AddPrimitiveData ( DataElement const & de )`

27.179.5.2 `static MPTYPE gdcm::MeshPrimitive::GetMPTYPE ( const char * type )` `[static]`

27.179.5.3 `static const char* gdcm::MeshPrimitive::GetMPTYPEString ( const MPTYPE type )` `[static]`

27.179.5.4 `unsigned int gdcm::MeshPrimitive::GetNumberOfPrimitivesData ( ) const`

27.179.5.5 `const DataElement& gdcm::MeshPrimitive::GetPrimitiveData ( ) const`

27.179.5.6 `DataElement& gdcm::MeshPrimitive::GetPrimitiveData ( )`

27.179.5.7 `const DataElement& gdcm::MeshPrimitive::GetPrimitiveData ( const unsigned int idx ) const`

27.179.5.8 `DataElement& gdcm::MeshPrimitive::GetPrimitiveData ( const unsigned int idx )`

27.179.5.9 `const PrimitivesData& gdcm::MeshPrimitive::GetPrimitivesData ( ) const`

27.179.5.10 `PrimitivesData& gdcm::MeshPrimitive::GetPrimitivesData ( )`

27.179.5.11 `MPTYPE gdcm::MeshPrimitive::GetPrimitiveType ( ) const`

27.179.5.12 `void gdcm::MeshPrimitive::SetPrimitiveData ( DataElement const & de )`

27.179.5.13 `void gdcm::MeshPrimitive::SetPrimitiveData ( const unsigned int idx, DataElement const & de )`

27.179.5.14 `void gdcm::MeshPrimitive::SetPrimitivesData ( PrimitivesData const & DEs )`

27.179.5.15 `void gdcm::MeshPrimitive::SetPrimitiveType ( const MPTYPE type )`

### 27.179.6 Member Data Documentation

27.179.6.1 `PrimitivesData gdcm::MeshPrimitive::PrimitiveData` `[protected]`

27.179.6.2 `MPTYPE gdcm::MeshPrimitive::PrimitiveType` `[protected]`


The documentation for this class was generated from the following file:

- [gdcmMeshPrimitive.h](#)


## 27.180 gdcm::ModifiedEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::ModifiedEvent`:

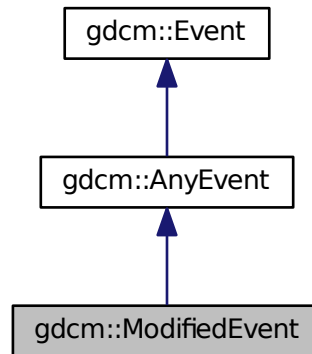

Collaboration diagram for `gdcm::ModifiedEvent`:

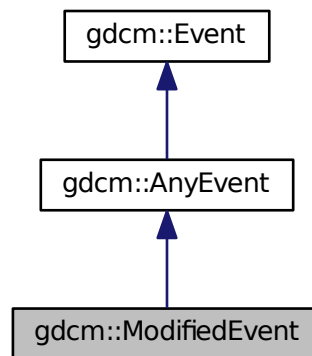

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.181 gdcm::Module Class Reference

Class for representing a [Module](#).

```
#include <gdcmModule.h>
```

### Public Types

- typedef std::vector< std::string > [ArrayIncludeMacrosType](#)
- typedef std::map< [Tag](#), [ModuleEntry](#) > [MapModuleEntry](#)

### Public Member Functions

- [Module](#) ()
- void [AddMacro](#) (const char \*include)
- void [AddModuleEntry](#) (const [Tag](#) &tag, const [ModuleEntry](#) &module)  
*Will add a [ModuleEntry](#) directly at root-level. See [Macro](#) for nested-included level.*
- void [Clear](#) ()
- bool [FindModuleEntryInMacros](#) ([Macros](#) const &macros, const [Tag](#) &tag) const
- const [ModuleEntry](#) & [GetModuleEntryInMacros](#) ([Macros](#) const &macros, const [Tag](#) &tag) const
- const char \* [GetName](#) () const
- void [SetName](#) (const char \*name)
- bool [Verify](#) (const [DataSet](#) &ds, [Usage](#) const &usage) const

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Module](#) &\_val)

### 27.181.1 Detailed Description

Class for representing a [Module](#).

#### Note

[Module](#): A set of Attributes within an Information Entity or Normalized [IOD](#) which are logically related to each other.

#### See also

[Macro](#)

#### Examples:

[TraverseModules.cxx](#).

### 27.181.2 Member Typedef Documentation

27.181.2.1 typedef std::vector<std::string> gdcm::Module::ArrayIncludeMacrosType

27.181.2.2 typedef std::map<Tag, ModuleEntry> gdcm::Module::MapModuleEntry

### 27.181.3 Constructor & Destructor Documentation

27.181.3.1 `gdcmmodule::Module ( ) [inline]`

### 27.181.4 Member Function Documentation

27.181.4.1 `void gdcmmodule::AddMacro ( const char * include ) [inline]`

27.181.4.2 `void gdcmmodule::AddModuleEntry ( const Tag & tag, const ModuleEntry & module ) [inline]`

Will add a [ModuleEntry](#) directly at root-level. See [Macro](#) for nested-included level.

27.181.4.3 `void gdcmmodule::Clear ( ) [inline]`

27.181.4.4 `bool gdcmmodule::FindModuleEntryInMacros ( Macros const & macros, const Tag & tag ) const`

Find or Get a [ModuleEntry](#). [ModuleEntry](#) are either search are root-level or within nested-macro included in module.

Examples:

[TraverseModules.cxx](#).

27.181.4.5 `const ModuleEntry& gdcmmodule::GetModuleEntryInMacros ( Macros const & macros, const Tag & tag ) const`

Examples:

[TraverseModules.cxx](#).

27.181.4.6 `const char* gdcmmodule::GetName ( ) const [inline]`

27.181.4.7 `void gdcmmodule::SetName ( const char * name ) [inline]`

27.181.4.8 `bool gdcmmodule::Verify ( const DataSet & ds, Usage const & usage ) const`

### 27.181.5 Friends And Related Function Documentation

27.181.5.1 `std::ostream& operator<< ( std::ostream & _os, const Module & _val ) [friend]`

The documentation for this class was generated from the following file:


- [gdcmmodule.h](#)

## 27.182 gdcmmodule::ModuleEntry Class Reference


Class for representing a [ModuleEntry](#).

```
#include <gdcmmoduleentry.h>
```

Inheritance diagram for gdcm::ModuleEntry:

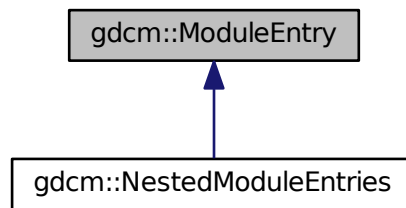

Collaboration diagram for gdcm::ModuleEntry:

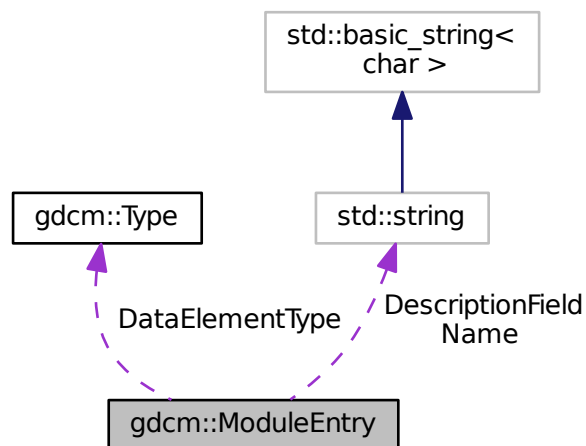

## Public Types

- typedef std::string [Description](#)

## Public Member Functions

- [ModuleEntry](#) (const char \*name="", const char \*type="3", const char \*description="")
- virtual [~ModuleEntry](#) ()
- const [Description](#) & [GetDescription](#) () const
- const char \* [GetName](#) () const

- const [Type](#) & [GetType](#) () const
- void [SetDescription](#) (const char \*d)
- void [SetName](#) (const char \*name)
- void [SetType](#) (const [Type](#) &type)

### Protected Attributes

- [Type](#) [DataElementType](#)
- [Description](#) [DescriptionField](#)
- std::string [Name](#)

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [ModuleEntry](#) &\_val)

## 27.182.1 Detailed Description

Class for representing a [ModuleEntry](#).

#### Note

bla

#### See also

[DictEntry](#)

#### Examples:

[TraverseModules.cxx](#).

## 27.182.2 Member Typedef Documentation

27.182.2.1 `typedef std::string gdcmm::ModuleEntry::Description`

## 27.182.3 Constructor & Destructor Documentation

27.182.3.1 `gdcmm::ModuleEntry::ModuleEntry ( const char * name = " ", const char * type = "3", const char * description = " " )`  
`[inline]`

References `gdcmm::Type::GetTypeType()`.

27.182.3.2 `virtual gdcmm::ModuleEntry::~~ModuleEntry ( )` `[inline]`, `[virtual]`

## 27.182.4 Member Function Documentation

27.182.4.1 `const Description& gdcmm::ModuleEntry::GetDescription ( ) const` `[inline]`

27.182.4.2 `const char* gdcmm::ModuleEntry::GetName ( ) const` `[inline]`



27.182.4.3 `const Type& gdcm::ModuleEntry::GetType ( ) const` `[inline]`

Examples:

[TraverseModules.cxx](#).

27.182.4.4 `void gdcm::ModuleEntry::SetDescription ( const char * d )` `[inline]`

27.182.4.5 `void gdcm::ModuleEntry::SetName ( const char * name )` `[inline]`

27.182.4.6 `void gdcm::ModuleEntry::SetType ( const Type & type )` `[inline]`

## 27.182.5 Friends And Related Function Documentation

27.182.5.1 `std::ostream& operator<< ( std::ostream & _os, const ModuleEntry & _val )` `[friend]`

## 27.182.6 Member Data Documentation

27.182.6.1 `Type gdcm::ModuleEntry::DataElementType` `[protected]`

Referenced by `gdcm::operator<<()`.

27.182.6.2 `Description gdcm::ModuleEntry::DescriptionField` `[protected]`

Referenced by `gdcm::operator<<()`.

27.182.6.3 `std::string gdcm::ModuleEntry::Name` `[protected]`

Referenced by `gdcm::operator<<()`.

The documentation for this class was generated from the following file:

- [gdcmModuleEntry.h](#)

## 27.183 gdcm::Modules Class Reference

Class for representing a [Modules](#).

```
#include <gdcmModules.h>
```

### Public Types

- `typedef std::map< std::string, Module > ModuleMapType`

### Public Member Functions

- [Modules](#) ()
- void [AddModule](#) (const char \*ref, const [Module](#) &module)

- void [Clear](#) ()
- const [Module](#) & [GetModule](#) (const char \*name) const
- bool [IsEmpty](#) () const

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Modules](#) &\_val)

### 27.183.1 Detailed Description

Class for representing a [Modules](#).

#### Note

bla

#### See also

[Module](#)

#### Examples:

[TraverseModules.cxx](#).

### 27.183.2 Member Typedef Documentation

27.183.2.1 `typedef std::map<std::string, Module> gdcm::Modules::ModuleMapType`

### 27.183.3 Constructor & Destructor Documentation

27.183.3.1 `gdcm::Modules::Modules ( )` `[inline]`

### 27.183.4 Member Function Documentation

27.183.4.1 `void gdcm::Modules::AddModule ( const char * ref, const Module & module )` `[inline]`

27.183.4.2 `void gdcm::Modules::Clear ( )` `[inline]`

27.183.4.3 `const Module& gdcm::Modules::GetModule ( const char * name ) const` `[inline]`

#### Examples:

[TraverseModules.cxx](#).

27.183.4.4 `bool gdcm::Modules::IsEmpty ( ) const` `[inline]`

### 27.183.5 Friends And Related Function Documentation

27.183.5.1 `std::ostream& operator<< ( std::ostream &_os, const Modules &_val )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmModules.h](#)

## 27.184 gdcm::MovePatientRootQuery Class Reference


[MovePatientRootQuery](#) contains: the class which will produce a dataset for c-move with patient root.

```
#include <gdcmMovePatientRootQuery.h>
```

Inheritance diagram for gdcm::MovePatientRootQuery:

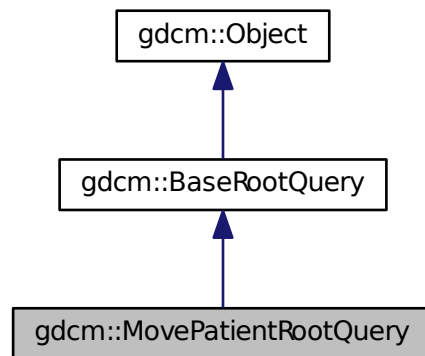

Collaboration diagram for gdcm::MovePatientRootQuery:

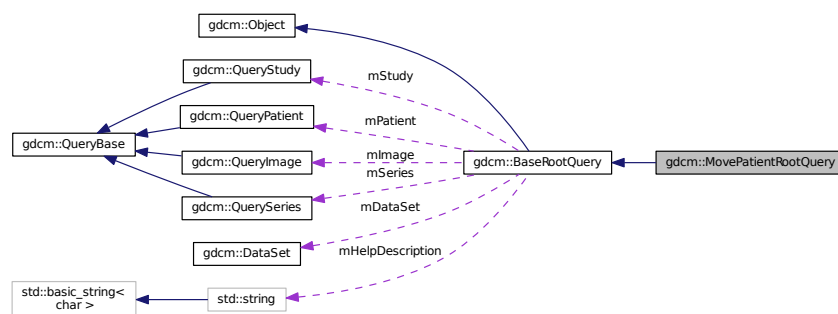

### Public Member Functions

- [MovePatientRootQuery](#) ()
- `UIDs::TSName GetAbstractSyntaxUID` () const
- `std::vector< Tag > GetTagListByLevel` (const [EQueryLevel](#) &inQueryLevel)
- void [InitializeDataSet](#) (const [EQueryLevel](#) &inQueryLevel)

- bool [ValidateQuery](#) (bool inStrict=true) const

## Friends

- class [QueryFactory](#)

## Additional Inherited Members

### 27.184.1 Detailed Description

[MovePatientRootQuery](#) contains: the class which will produce a dataset for c-move with patient root.

### 27.184.2 Constructor & Destructor Documentation

27.184.2.1 `gdcm::MovePatientRootQuery::MovePatientRootQuery ( )`

### 27.184.3 Member Function Documentation

27.184.3.1 `UIDs::TSName gdcm::MovePatientRootQuery::GetAbstractSyntaxUID ( ) const` [virtual]

Implements [gdcm::BaseRootQuery](#).

27.184.3.2 `std::vector<Tag> gdcm::MovePatientRootQuery::GetTagListByLevel ( const EQueryLevel & inQueryLevel )`  
[virtual]

this function will return all tags at a given query level, so that they maybe selected for searching. The boolean forFind is true if the query is a find query, or false for a move query.

Implements [gdcm::BaseRootQuery](#).

27.184.3.3 `void gdcm::MovePatientRootQuery::InitializeDataSet ( const EQueryLevel & inQueryLevel )` [virtual]

this function sets tag 8,52 to the appropriate value based on query level also fills in the right unique tags, as per the standard's requirements should allow for connection with dcm4chee

Implements [gdcm::BaseRootQuery](#).

27.184.3.4 `bool gdcm::MovePatientRootQuery::ValidateQuery ( bool inStrict=true ) const` [virtual]

have to be able to ensure that 0x8,0x52 is set (which will be true if InitializeDataSet is called...) that the level is appropriate (ie, not setting PATIENT for a study query that the tags in the query match the right level (either required, unique, optional) by default, this function checks to see if the query is for finding, which is more permissive than for moving. For moving, only the unique tags are allowed. 10 Jan 2011: adding in the 'strict' mode. according to the standard (at least, how I've read it), only tags for a particular level should be allowed in a particular query (ie, just series level tags in a series level query). However, it seems that dcm4chee doesn't share that interpretation. So, if 'inStrict' is false, then tags from the current level and all higher levels are now considered valid. So, if you're doing a non-strict series-level query, tags from the patient and study level can be passed along as well.

Implements [gdcm::BaseRootQuery](#).

## 27.184.4 Friends And Related Function Documentation

### 27.184.4.1 friend class QueryFactory [friend]

The documentation for this class was generated from the following file:


- [gdcmmovePatientRootQuery.h](#)

## 27.185 gdcmmoveStudyRootQuery Class Reference


[MoveStudyRootQuery](#) contains: the class which will produce a dataset for C-MOVE with study root.

```
#include <gdcmmoveStudyRootQuery.h>
```

Inheritance diagram for gdcmmoveStudyRootQuery:

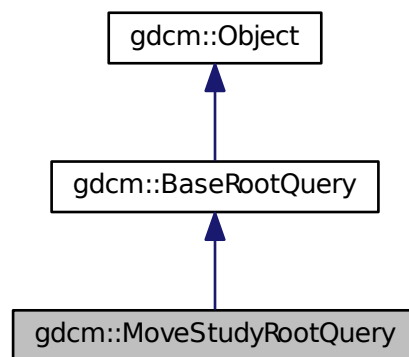

Collaboration diagram for gdcmmoveStudyRootQuery:

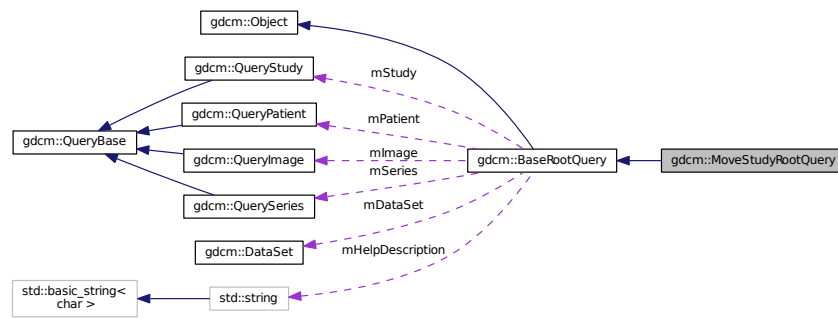

## Public Member Functions

- [MoveStudyRootQuery](#) ()
- [UIDs::TSName GetAbstractSyntaxUID](#) () const
- `std::vector< Tag > GetTagListByLevel (const EQueryLevel &inQueryLevel)`
- `void InitializeDataSet (const EQueryLevel &inQueryLevel)`
- `bool ValidateQuery (bool inStrict=true) const`

## Friends

- class [QueryFactory](#)

## Additional Inherited Members

### 27.185.1 Detailed Description

[MoveStudyRootQuery](#) contains: the class which will produce a dataset for C-MOVE with study root.

### 27.185.2 Constructor & Destructor Documentation

27.185.2.1 `gdcm::MoveStudyRootQuery::MoveStudyRootQuery ( )`

### 27.185.3 Member Function Documentation

27.185.3.1 `UIDs::TSName gdcm::MoveStudyRootQuery::GetAbstractSyntaxUID ( ) const` `[virtual]`

Implements [gdcm::BaseRootQuery](#).

27.185.3.2 `std::vector<Tag> gdcm::MoveStudyRootQuery::GetTagListByLevel ( const EQueryLevel & inQueryLevel )`  
`[virtual]`

this function will return all tags at a given query level, so that they maybe selected for searching. The boolean forFind is true if the query is a find query, or false for a move query.

Implements [gdcm::BaseRootQuery](#).

27.185.3.3 `void gdcm::MoveStudyRootQuery::InitializeDataSet ( const EQueryLevel & inQueryLevel )` `[virtual]`

this function sets tag 8,52 to the appropriate value based on query level also fills in the right unique tags, as per the standard's requirements should allow for connection with dcmTk

Implements [gdcm::BaseRootQuery](#).

27.185.3.4 `bool gdcm::MoveStudyRootQuery::ValidateQuery ( bool inStrict = true ) const` `[virtual]`

have to be able to ensure that 0x8,0x52 is set (which will be true if InitializeDataSet is called...) that the level is appropriate (ie, not setting PATIENT for a study query that the tags in the query match the right level (either required, unique, optional) by default, this function checks to see if the query is for finding, which is more permissive than for moving. For moving, only the unique tags are allowed. 10 Jan 2011: adding in the 'strict' mode. according to the

standard (at least, how I've read it), only tags for a particular level should be allowed in a particular query (ie, just series level tags in a series level query). However, it seems that dcm4chee doesn't share that interpretation. So, if 'inStrict' is false, then tags from the current level and all higher levels are now considered valid. So, if you're doing a non-strict series-level query, tags from the patient and study level can be passed along as well.

Implements [gdcm::BaseRootQuery](#).

## 27.185.4 Friends And Related Function Documentation

### 27.185.4.1 friend class QueryFactory [friend]

The documentation for this class was generated from the following file:


- [gdcmMoveStudyRootQuery.h](#)

## 27.186 gdcm::NestedModuleEntries Class Reference


Class for representing a [NestedModuleEntries](#).

```
#include <gdcmNestedModuleEntries.h>
```

Inheritance diagram for gdcm::NestedModuleEntries:

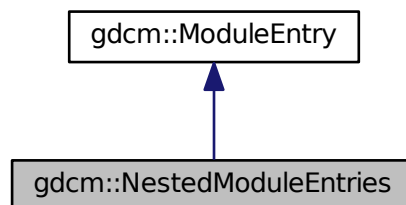

Collaboration diagram for `gdcm::NestedModuleEntries`:

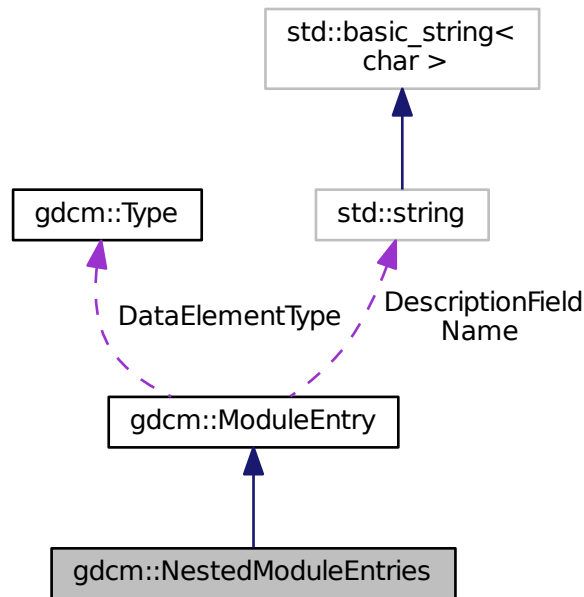

## Public Types

- typedef `std::vector< ModuleEntry >::size_type` [SizeType](#)

## Public Member Functions

- [NestedModuleEntries](#) (const char \*name="", const char \*type="3", const char \*description="")
- void [AddModuleEntry](#) (const [ModuleEntry](#) &me)
- const [ModuleEntry](#) & [GetModuleEntry](#) ([SizeType](#) idx) const
- [ModuleEntry](#) & [GetModuleEntry](#) ([SizeType](#) idx)
- [SizeType](#) [GetNumberOfModuleEntries](#) ()

## Friends

- `std::ostream & operator<< (std::ostream &_os, const NestedModuleEntries &_val)`

## Additional Inherited Members

### 27.186.1 Detailed Description

Class for representing a [NestedModuleEntries](#).



## Note

bla

## See also

[ModuleEntry](#)

## 27.186.2 Member Typedef Documentation

27.186.2.1 `typedef std::vector<ModuleEntry>::size_type gdcm::NestedModuleEntries::SizeType`

## 27.186.3 Constructor &amp; Destructor Documentation

27.186.3.1 `gdcm::NestedModuleEntries::NestedModuleEntries ( const char * name = " ", const char * type = "3", const char * description = " " )` `[inline]`

## 27.186.4 Member Function Documentation

27.186.4.1 `void gdcm::NestedModuleEntries::AddModuleEntry ( const ModuleEntry & me )` `[inline]`

27.186.4.2 `const ModuleEntry& gdcm::NestedModuleEntries::GetModuleEntry ( SizeType idx ) const` `[inline]`

27.186.4.3 `ModuleEntry& gdcm::NestedModuleEntries::GetModuleEntry ( SizeType idx )` `[inline]`

27.186.4.4 `SizeType gdcm::NestedModuleEntries::GetNumberOfModuleEntries ( )` `[inline]`

## 27.186.5 Friends And Related Function Documentation

27.186.5.1 `std::ostream& operator<< ( std::ostream & _os, const NestedModuleEntries & _val )` `[friend]`


The documentation for this class was generated from the following file:

- [gdcmNestedModuleEntries.h](#)


## 27.187 gdcm::NoEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::NoEvent`:

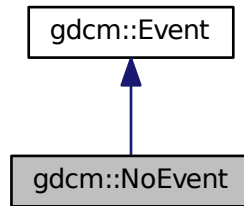

Collaboration diagram for `gdcm::NoEvent`:

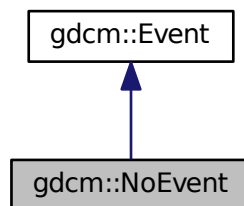

### Additional Inherited Members

#### 27.187.1 Detailed Description

Define some common GDCM events

The documentation for this class was generated from the following file:


- [gdcmEvent.h](#)

## 27.188 `gdcm::Object` Class Reference

[Object.](#)

```
#include <gdcmObject.h>
```

Inheritance diagram for gdcm::Object:

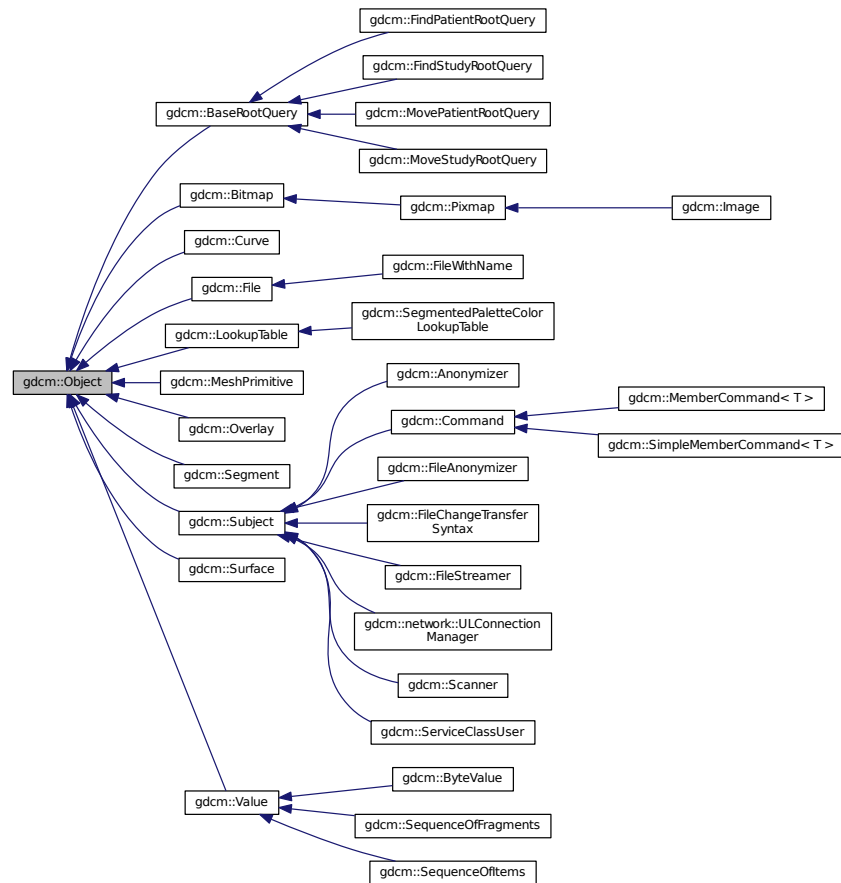

## Public Member Functions

- [Object](#) ()
- [Object](#) (const [Object](#) &)
- *Special requirement for copy/cstor, assignment operator.*
- virtual [~Object](#) ()
- void [operator=](#) (const [Object](#) &)
- virtual void [Print](#) (std::ostream &) const

## Protected Member Functions

- void [Register](#) ()
- void [UnRegister](#) ()

## Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [Object](#) &obj)

- `template<class ObjectType >`  
class [SmartPointer](#)

### 27.188.1 Detailed Description

[Object](#).

Note

main superclass for object that want to use [SmartPointer](#) invasive ref counting system

See also

[SmartPointer](#)

### 27.188.2 Constructor & Destructor Documentation

27.188.2.1 `gdcM::Object::Object ( )` `[inline]`

27.188.2.2 `virtual gdcM::Object::~~Object ( )` `[inline], [virtual]`

27.188.2.3 `gdcM::Object::Object ( const Object & )` `[inline]`

Special requirement for copy/cstor, assignment operator.

### 27.188.3 Member Function Documentation

27.188.3.1 `void gdcM::Object::operator= ( const Object & )` `[inline]`

27.188.3.2 `virtual void gdcM::Object::Print ( std::ostream & ) const` `[inline], [virtual]`

Reimplemented in [gdcM::SequenceOfFragments](#), [gdcM::ByteValue](#), [gdcM::SequenceOfItems](#), [gdcM::BaseRootQuery](#), [gdcM::Scanner](#), [gdcM::Image](#), [gdcM::Curve](#), [gdcM::Overlay](#), [gdcM::Bitmap](#), [gdcM::LookupTable](#), [gdcM::Pixmap](#), and [gdcM::SegmentedPaletteColorLookupTable](#).

Examples:

[ReadAndDumpDICOMDIR.cxx](#).

Referenced by `gdcM::operator<<()`.

27.188.3.3 `void gdcM::Object::Register ( )` `[inline], [protected]`

27.188.3.4 `void gdcM::Object::UnRegister ( )` `[inline], [protected]`

### 27.188.4 Friends And Related Function Documentation

27.188.4.1 `std::ostream& operator<< ( std::ostream & os, const Object & obj )` `[friend]`

27.188.4.2 `template<class ObjectType > friend class SmartPointer` [`friend`]


The documentation for this class was generated from the following file:

- [gdcMObject.h](#)


## 27.189 gdcM::OpenSSLCryptoFactory Class Reference

```
#include <gdcMOpenSSLCryptoFactory.h>
```

Inheritance diagram for gdcM::OpenSSLCryptoFactory:

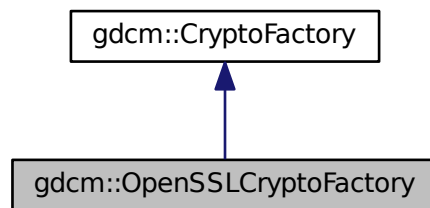

Collaboration diagram for gdcM::OpenSSLCryptoFactory:

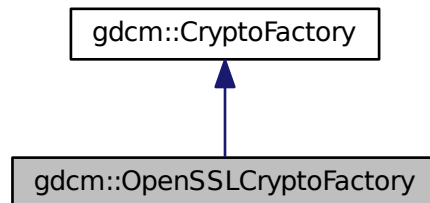

### Public Member Functions

- [OpenSSLCryptoFactory](#) ([CryptoLib](#) id)
- [CryptographicMessageSyntax](#) \* [CreateCMSProvider](#) ()

### Protected Member Functions

- void [InitOpenSSL](#) ()

## Additional Inherited Members

### 27.189.1 Constructor & Destructor Documentation

27.189.1.1 `gdcmm::OpenSSLCryptoFactory::OpenSSLCryptoFactory ( CryptoLib id )` `[inline]`

References `gdcmDebugMacro`.

### 27.189.2 Member Function Documentation

27.189.2.1 `CryptographicMessageSyntax*` `gdcmm::OpenSSLCryptoFactory::CreateCMSProvider ( )` `[inline]`,  
`[virtual]`

Implements [gdcmm::CryptoFactory](#).

27.189.2.2 `void` `gdcmm::OpenSSLCryptoFactory::InitOpenSSL ( )` `[protected]`


The documentation for this class was generated from the following file:

- [gdcmmOpenSSLCryptoFactory.h](#)


## 27.190 gdcmm::OpenSSLCryptographicMessageSyntax Class Reference

```
#include <gdcmmOpenSSLCryptographicMessageSyntax.h>
```

Inheritance diagram for `gdcmm::OpenSSLCryptographicMessageSyntax`:

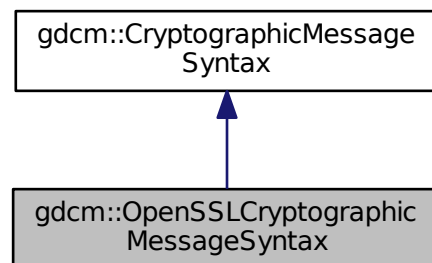

Collaboration diagram for gdcmm::OpenSSLCryptographicMessageSyntax:

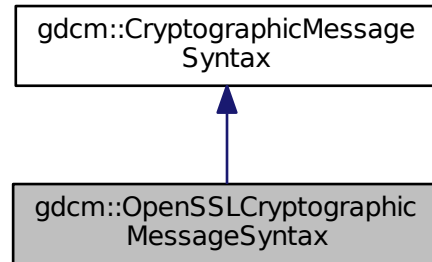

## Public Member Functions

- [OpenSSLCryptographicMessageSyntax \(\)](#)
- [~OpenSSLCryptographicMessageSyntax \(\)](#)
- bool [Decrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*decrypt content from a PKCS#7 envelopedData structure*
- bool [Encrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*create a CMS envelopedData structure*
- [CipherTypes GetCipherType](#) () const
- bool [ParseCertificateFile](#) (const char \*filename)
- bool [ParseKeyFile](#) (const char \*filename)
- void [SetCipherType](#) ([CipherTypes](#) type)
- bool [SetPassword](#) (const char \*pass, size\_t passLen)

## Additional Inherited Members

### 27.190.1 Constructor & Destructor Documentation

27.190.1.1 `gdcmm::OpenSSLCryptographicMessageSyntax::OpenSSLCryptographicMessageSyntax ( )`

27.190.1.2 `gdcmm::OpenSSLCryptographicMessageSyntax::~~OpenSSLCryptographicMessageSyntax ( )`

### 27.190.2 Member Function Documentation

27.190.2.1 `bool gdcmm::OpenSSLCryptographicMessageSyntax::Decrypt ( char * output, size_t & outlen, const char * array, size_t len ) const` [virtual]

decrypt content from a PKCS#7 envelopedData structure

Implements [gdcmm::CryptographicMessageSyntax](#).

27.190.2.2 `bool gdcM::OpenSSLCryptographicMessageSyntax::Encrypt ( char * output, size_t & outlen, const char * array, size_t len ) const` `[virtual]`

create a CMS envelopedData structure

Implements [gdcM::CryptographicMessageSyntax](#).

27.190.2.3 `CipherTypes gdcM::OpenSSLCryptographicMessageSyntax::GetCipherType ( ) const` `[virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.190.2.4 `bool gdcM::OpenSSLCryptographicMessageSyntax::ParseCertificateFile ( const char * filename )` `[virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.190.2.5 `bool gdcM::OpenSSLCryptographicMessageSyntax::ParseKeyFile ( const char * filename )` `[virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

27.190.2.6 `void gdcM::OpenSSLCryptographicMessageSyntax::SetCipherType ( CipherTypes type )` `[virtual]`

Set Cipher [Type](#). Default is: AES256\_CIPHER

Implements [gdcM::CryptographicMessageSyntax](#).

27.190.2.7 `bool gdcM::OpenSSLCryptographicMessageSyntax::SetPassword ( const char * pass, size_t passLen )` `[virtual]`

Implements [gdcM::CryptographicMessageSyntax](#).

The documentation for this class was generated from the following file:


- [gdcMOpenSSLCryptographicMessageSyntax.h](#)

## 27.191 gdcM::OpenSSLP7CryptoFactory Class Reference


```
#include <gdcMOpenSSLP7CryptoFactory.h>
```



Inheritance diagram for gdcmm::OpenSSLP7CryptoFactory:

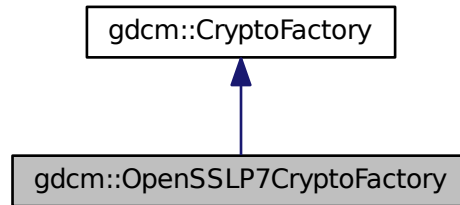

Collaboration diagram for gdcmm::OpenSSLP7CryptoFactory:

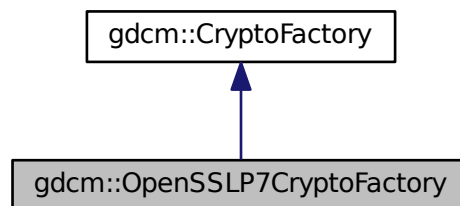

## Public Member Functions

- [OpenSSLP7CryptoFactory](#) ([CryptoLib](#) id)
- [CryptographicMessageSyntax](#) \* [CreateCMSProvider](#) ()

## Additional Inherited Members

### 27.191.1 Constructor & Destructor Documentation

27.191.1.1 `gdcmm::OpenSSLP7CryptoFactory::OpenSSLP7CryptoFactory ( CryptoLib id )` `[inline]`

References `gdcmmDebugMacro`.

### 27.191.2 Member Function Documentation

27.191.2.1 **CryptographicMessageSyntax\*** `gdcM::OpenSSLP7CryptoFactory::CreateCMSProvider ( )` `[inline]`,  
`[virtual]`

Implements [gdcM::CryptoFactory](#).

The documentation for this class was generated from the following file:


- [gdcMOpenSSLP7CryptoFactory.h](#)

## 27.192 gdcM::OpenSSLP7CryptographicMessageSyntax Class Reference


Class for [CryptographicMessageSyntax](#) encryption. This is just a simple wrapper around openssl PKCS7\_encrypt functionalities.

```
#include <gdcMOpenSSLP7CryptographicMessageSyntax.h>
```

Inheritance diagram for `gdcM::OpenSSLP7CryptographicMessageSyntax`:

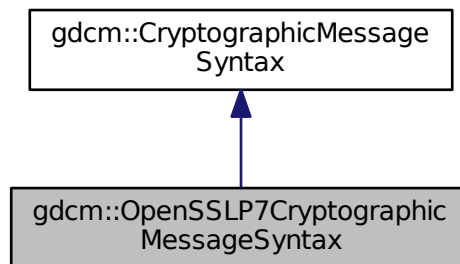

Collaboration diagram for `gdcM::OpenSSLP7CryptographicMessageSyntax`:

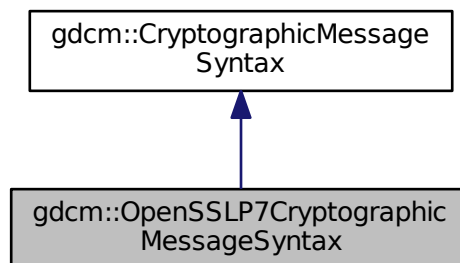

## Public Member Functions

- [OpenSSLP7CryptographicMessageSyntax](#) ()
- [~OpenSSLP7CryptographicMessageSyntax](#) ()
- bool [Decrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*decrypt content from a PKCS#7 envelopedData structure*
- bool [Encrypt](#) (char \*output, size\_t &outlen, const char \*array, size\_t len) const  
*create a PKCS#7 envelopedData structure*
- [CipherTypes](#) [GetCipherType](#) () const
- bool [ParseCertificateFile](#) (const char \*filename)
- bool [ParseKeyFile](#) (const char \*filename)
- void [SetCipherType](#) ([CipherTypes](#) type)
- bool [SetPassword](#) (const char \*, size\_t)

## Additional Inherited Members

### 27.192.1 Detailed Description

Class for [CryptographicMessageSyntax](#) encryption. This is just a simple wrapper around openssl PKCS7\_encrypt functionalities.

See online documentation [http://www.openssl.org/docs/crypto/PKCS7\\_encrypt.html](http://www.openssl.org/docs/crypto/PKCS7_encrypt.html)

### 27.192.2 Constructor & Destructor Documentation

27.192.2.1 gdcmm::OpenSSLP7CryptographicMessageSyntax::OpenSSLP7CryptographicMessageSyntax ( )

27.192.2.2 gdcmm::OpenSSLP7CryptographicMessageSyntax::~~OpenSSLP7CryptographicMessageSyntax ( )

### 27.192.3 Member Function Documentation

27.192.3.1 bool gdcmm::OpenSSLP7CryptographicMessageSyntax::Decrypt ( char \* output, size\_t & outlen, const char \* array, size\_t len ) const [virtual]

decrypt content from a PKCS#7 envelopedData structure

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.2 bool gdcmm::OpenSSLP7CryptographicMessageSyntax::Encrypt ( char \* output, size\_t & outlen, const char \* array, size\_t len ) const [virtual]

create a PKCS#7 envelopedData structure

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.3 [CipherTypes](#) gdcmm::OpenSSLP7CryptographicMessageSyntax::GetCipherType ( ) const [virtual]

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.4 `bool gdcmm::OpenSSL7CryptographicMessageSyntax::ParseCertificateFile ( const char * filename )` [virtual]

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.5 `bool gdcmm::OpenSSL7CryptographicMessageSyntax::ParseKeyFile ( const char * filename )` [virtual]

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.6 `void gdcmm::OpenSSL7CryptographicMessageSyntax::SetCipherType ( CipherTypes type )` [virtual]

Set Cipher [Type](#). Default is: AES256\_CIPHER

Implements [gdcmm::CryptographicMessageSyntax](#).

27.192.3.7 `bool gdcmm::OpenSSL7CryptographicMessageSyntax::SetPassword ( const char *, size_t )` [inline],  
[virtual]

Implements [gdcmm::CryptographicMessageSyntax](#).

References `gdcmmWarningMacro`.

The documentation for this class was generated from the following file:

- [gdcmmOpenSSL7CryptographicMessageSyntax.h](#)

## 27.193 gdcmm::Orientation Class Reference

class to handle [Orientation](#)

```
#include <gdcmmOrientation.h>
```

### Public Types

- enum [OrientationType](#) {  
    [UNKNOWN](#),  
    [AXIAL](#),  
    [CORONAL](#),  
    [SAGITTAL](#),  
    [OBLIQUE](#) }

### Public Member Functions

- [Orientation](#) ()
- [~Orientation](#) ()
- void [Print](#) (std::ostream &) const

*Print.*

## Static Public Member Functions

- static const char \* [GetLabel](#) ([OrientationType](#) type)  
*Return the label of an [Orientation](#).*
- static double [GetObliquityThresholdCosineValue](#) ()
- static [OrientationType](#) [GetType](#) (const double dircos[6])
- static void [SetObliquityThresholdCosineValue](#) (double val)  
*ObliquityThresholdCosineValue stuff.*

## Static Protected Member Functions

- static char [GetMajorAxisFromPatientRelativeDirectionCosine](#) (double x, double y, double z)

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Orientation](#) &o)

### 27.193.1 Detailed Description

class to handle [Orientation](#)

### 27.193.2 Member Enumeration Documentation

#### 27.193.2.1 enum gdcm::Orientation::OrientationType

Enumerator

**UNKNOWN**

**AXIAL**

**CORONAL**

**SAGITTAL**

**OBLIQUE**

### 27.193.3 Constructor & Destructor Documentation

#### 27.193.3.1 gdcm::Orientation::Orientation ( )

#### 27.193.3.2 gdcm::Orientation::~~Orientation ( )

### 27.193.4 Member Function Documentation

#### 27.193.4.1 static const char\* gdcm::Orientation::GetLabel ( [OrientationType](#) type ) [static]

Return the label of an [Orientation](#).

27.193.4.2 `static char gdcM::Orientation::GetMajorAxisFromPatientRelativeDirectionCosine ( double x, double y, double z )`  
`[static], [protected]`

27.193.4.3 `static double gdcM::Orientation::GetObliquityThresholdCosineValue ( )` `[static]`

27.193.4.4 `static OrientationType gdcM::Orientation::GetType ( const double dircos[6] )` `[static]`

Return the type of orientation from a direction cosines Input is an array of 6 double

27.193.4.5 `void gdcM::Orientation::Print ( std::ostream & ) const`

Print.

Referenced by `gdcM::operator<<()`.

27.193.4.6 `static void gdcM::Orientation::SetObliquityThresholdCosineValue ( double val )` `[static]`

ObliquityThresholdCosineValue stuff.

## 27.193.5 Friends And Related Function Documentation

27.193.5.1 `std::ostream& operator<< ( std::ostream &_os, const Orientation &o )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcMOrientation.h](#)

## 27.194 gdcM::Overlay Class Reference


[Overlay](#) class.

```
#include <gdcMOverlay.h>
```

Inheritance diagram for `gdcM::Overlay`:

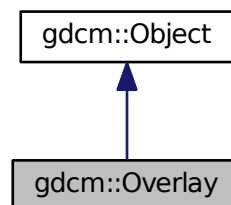

Collaboration diagram for gdcm::Overlay:

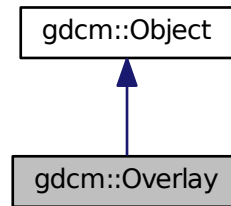

## Public Types

- enum `OverlayType` {  
     `Invalid` = 0,  
     `Graphics` = 1,  
     `ROI` = 2 }

## Public Member Functions

- `Overlay ()`
- `Overlay (Overlay const &ov)`
- `~Overlay ()`
- void `Decompress (std::ostream &os) const`  
     *Decode the internal OverlayData (packed bits) into unpacked representation.*
- unsigned short `GetBitPosition () const`  
     *return bit position*
- unsigned short `GetBitsAllocated () const`  
     *return bits allocated*
- unsigned short `GetColumns () const`  
     *get columns*
- const char \* `GetDescription () const`  
     *get description*
- unsigned short `GetGroup () const`  
     *Get Group number.*
- const signed short \* `GetOrigin () const`  
     *get origin*
- const `ByteValue` & `GetOverlayData () const`
- unsigned short `GetRows () const`  
     *get rows*
- const char \* `GetType () const`  
     *get type*
- `OverlayType` `GetTypeAsEnum () const`

- bool [GetUnpackBuffer](#) (char \*buffer, size\_t len) const
- size\_t [GetUnpackBufferLength](#) () const
- bool [GrabOverlayFromPixelData](#) ([DataSet](#) const &ds)
- bool [IsEmpty](#) () const  
*Return whether or not the [Overlay](#) is empty:*
- bool [IsInPixelData](#) () const  
*return if the [Overlay](#) is stored in the pixel data or not*
- void [IsInPixelData](#) (bool b)  
*Set whether or not the OverlayData is in the Pixel Data:*
- bool [IsZero](#) () const  
*return true if all bits are set to 0*
- void [Print](#) (std::ostream &) const  
*Print.*
- void [SetBitPosition](#) (unsigned short bitposition)  
*set bit position*
- void [SetBitsAllocated](#) (unsigned short bitsallocated)  
*set bits allocated*
- void [SetColumns](#) (unsigned short columns)  
*set columns*
- void [SetDescription](#) (const char \*description)  
*set description*
- void [SetFrameOrigin](#) (unsigned short frameorigin)  
*set frame origin*
- void [SetGroup](#) (unsigned short group)  
*Set Group number.*
- void [SetNumberOfFrames](#) (unsigned int numberofframes)  
*set number of frames*
- void [SetOrigin](#) (const signed short origin[2])  
*set origin*
- void [SetOverlay](#) (const char \*array, size\_t length)  
*set overlay from byte array + length*
- void [SetRows](#) (unsigned short rows)  
*set rows*
- void [SetType](#) (const char \*type)  
*set type*
- void [Update](#) (const [DataElement](#) &de)  
*Update overlay from data element de:*

## Static Public Member Functions

- static const char \* [GetOverlayTypeAsString](#) ([OverlayType](#) ot)
- static [OverlayType](#) [GetOverlayTypeFromString](#) (const char \*)



## Additional Inherited Members

### 27.194.1 Detailed Description

[Overlay](#) class.

Note

see AreOverlaysInPixelData

**Todo** Is there actually any way to recognize an overlay ? On images with multiple overlay I do not see any way to differentiate them (other than the group tag).

Example:

### 27.194.2 Member Enumeration Documentation

#### 27.194.2.1 enum gdcm::Overlay::OverlayType

Enumerator

***Invalid***  
***Graphics***  
***ROI***

### 27.194.3 Constructor & Destructor Documentation

27.194.3.1 gdcm::Overlay::Overlay ( )

27.194.3.2 gdcm::Overlay::~~Overlay ( )

27.194.3.3 gdcm::Overlay::Overlay ( Overlay const & ov )

### 27.194.4 Member Function Documentation

27.194.4.1 void gdcm::Overlay::Decompress ( std::ostream & os ) const

Decode the internal OverlayData (packed bits) into unpacked representation.

27.194.4.2 unsigned short gdcm::Overlay::GetBitPosition ( ) const

return bit position

27.194.4.3 unsigned short gdcm::Overlay::GetBitsAllocated ( ) const

return bits allocated

27.194.4.4 unsigned short gdcm::Overlay::GetColumns ( ) const

get columns

27.194.4.5 `const char* gdcm::Overlay::GetDescription ( ) const`

get description

27.194.4.6 `unsigned short gdcm::Overlay::GetGroup ( ) const`

Get Group number.

27.194.4.7 `const signed short* gdcm::Overlay::GetOrigin ( ) const`

get origin

27.194.4.8 `const ByteValue& gdcm::Overlay::GetOverlayData ( ) const`

Return the [Overlay](#) Data as [ByteValue](#): Not thread safe

27.194.4.9 `static const char* gdcm::Overlay::GetOverlayTypeAsString ( OverlayType ot ) [static]`

27.194.4.10 `static OverlayType gdcm::Overlay::GetOverlayTypeFromString ( const char * ) [static]`

27.194.4.11 `unsigned short gdcm::Overlay::GetRows ( ) const`

get rows

27.194.4.12 `const char* gdcm::Overlay::GetType ( ) const`

get type

27.194.4.13 `OverlayType gdcm::Overlay::GetTypeAsEnum ( ) const`

27.194.4.14 `bool gdcm::Overlay::GetUnpackBuffer ( char * buffer, size_t len ) const`

Retrieve the unpack buffer for [Overlay](#). This is an error if the size if below [GetUnpackBufferLength\(\)](#)

27.194.4.15 `size_t gdcm::Overlay::GetUnpackBufferLength ( ) const`

Retrieve the size of the buffer needed to hold the [Overlay](#) as specified by Col & Row parameters

27.194.4.16 `bool gdcm::Overlay::GrabOverlayFromPixelData ( DataSet const & ds )`

27.194.4.17 `bool gdcm::Overlay::IsEmpty ( ) const`

Return whether or not the [Overlay](#) is empty:

27.194.4.18 `bool gdcm::Overlay::IsInPixelData ( ) const`

return if the [Overlay](#) is stored in the pixel data or not

27.194.4.19 void gdcmm::Overlay::IsInPixelData ( bool *b* )

Set whether or no the OverlayData is in the Pixel Data:

27.194.4.20 bool gdcmm::Overlay::IsZero ( ) const

return true if all bits are set to 0

27.194.4.21 void gdcmm::Overlay::Print ( std::ostream & ) const [virtual]

Print.

Reimplemented from [gdcmm::Object](#).

27.194.4.22 void gdcmm::Overlay::SetBitPosition ( unsigned short *bitposition* )

set bit position

27.194.4.23 void gdcmm::Overlay::SetBitsAllocated ( unsigned short *bitsallocated* )

set bits allocated

27.194.4.24 void gdcmm::Overlay::SetColumns ( unsigned short *columns* )

set columns

27.194.4.25 void gdcmm::Overlay::SetDescription ( const char \* *description* )

set description

27.194.4.26 void gdcmm::Overlay::SetFrameOrigin ( unsigned short *frameorigin* )

set frame origin

27.194.4.27 void gdcmm::Overlay::SetGroup ( unsigned short *group* )

Set Group number.

27.194.4.28 void gdcmm::Overlay::SetNumberOfFrames ( unsigned int *numberofframes* )

set number of frames

27.194.4.29 void gdcmm::Overlay::SetOrigin ( const signed short *origin*[2] )

set origin

27.194.4.30 void `gdcm::Overlay::SetOverlay` ( const char \* *array*, size\_t *length* )

set overlay from byte array + length

27.194.4.31 void `gdcm::Overlay::SetRows` ( unsigned short *rows* )

set rows

27.194.4.32 void `gdcm::Overlay::SetType` ( const char \* *type* )

set type

27.194.4.33 void `gdcm::Overlay::Update` ( const `DataElement` & *de* )

Update overlay from data element de:

The documentation for this class was generated from the following file:


- [gdcmOverlay.h](#)

## 27.195 gdcm::ParseException Class Reference


[ParseException](#) Standard exception handling object.

```
#include <gdcmParseException.h>
```

Inheritance diagram for `gdcm::ParseException`:

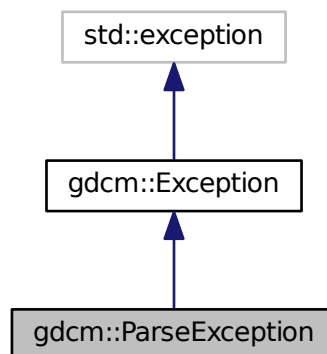

Collaboration diagram for gdcm::ParseException:

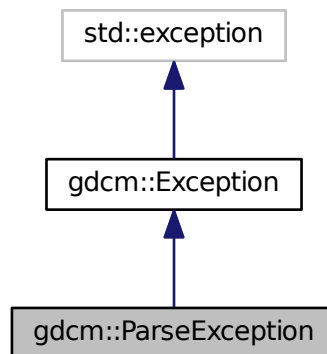

## Public Member Functions

- [ParseException](#) ()
- virtual [~ParseException](#) () throw ()
- const [DataElement](#) & [GetLastElement](#) () const
- [ParseException](#) & [operator=](#) (const [ParseException](#) &orig)
- void [SetLastElement](#) ([DataElement](#) &de)

### 27.195.1 Detailed Description

[ParseException](#) Standard exception handling object.

### 27.195.2 Constructor & Destructor Documentation

27.195.2.1 `gdcm::ParseException::ParseException ( )` [\[inline\]](#)

27.195.2.2 `virtual gdcm::ParseException::~~ParseException ( ) throw ()` [\[inline\]](#), [\[virtual\]](#)

### 27.195.3 Member Function Documentation

27.195.3.1 `const DataElement& gdcm::ParseException::GetLastElement ( ) const` [\[inline\]](#)

27.195.3.2 `ParseException& gdcm::ParseException::operator= ( const ParseException & orig )` [\[inline\]](#)

Assignment operator.

27.195.3.3 `void gdcm::ParseException::SetLastElement ( DataElement & de )` [\[inline\]](#)

Equivalence operator.

Referenced by `gdcm::Fragment::ReadBacktrack()`, and `gdcm::Fragment::ReadValue()`.

The documentation for this class was generated from the following file:

- [gdcmParseException.h](#)

## 27.196 gdcm::Parser Class Reference

[Parser](#) ala XML\_Parser from expat (SAX)

```
#include <gdcmParser.h>
```

### Public Types

- typedef void(\* [EndElementHandler](#)) (void \*userData, const [Tag](#) &name)
- enum [ErrorType](#) {  
[NoError](#),  
[NoMemoryError](#),  
[SyntaxError](#),  
[NoElementsError](#),  
[TagMismatchError](#),  
[DuplicateAttributeError](#),  
[JunkAfterDocElementError](#),  
[UndefinedEntityError](#),  
[UnexpectedStateError](#) }
- typedef void(\* [StartElementHandler](#)) (void \*userData, const [Tag](#) &tag, const char \*atts[])

### Public Member Functions

- [Parser](#) ()
- [~Parser](#) ()
- unsigned long [GetCurrentByteIndex](#) () const
- [ErrorType](#) [GetErrorCode](#) () const
- void \* [GetUserData](#) () const
- bool [Parse](#) (const char \*s, int len, bool isFinal)
- void [SetElementHandler](#) ([StartElementHandler](#) start, [EndElementHandler](#) end)
- void [SetUserData](#) (void \*userData)

### Static Public Member Functions

- static const char \* [GetErrorString](#) ([ErrorType](#) const &err)

### Protected Member Functions

- char \* [GetBuffer](#) (int len)
- bool [ParseBuffer](#) (int len, bool isFinal)
- [ErrorType](#) [Process](#) ()

### 27.196.1 Detailed Description

[Parser](#) ala XML\_Parser from expat (SAX)

Detailed description here

Note

Simple API for DICOM

### 27.196.2 Member Typedef Documentation

27.196.2.1 `typedef void(* gdcm::Parser::EndElementHandler) (void *userData, const Tag &name)`

27.196.2.2 `typedef void(* gdcm::Parser::StartElementHandler) (void *userData, const Tag &tag, const char *atts[ ])`

### 27.196.3 Member Enumeration Documentation

27.196.3.1 `enum gdcm::Parser::ErrorType`

Enumerator

***NoError***

***NoMemoryError***

***SyntaxError***

***NoElementsError***

***TagMismatchError***

***DuplicateAttributeError***

***JunkAfterDocElementError***

***UndefinedEntityError***

***UnexpectedStateError***

### 27.196.4 Constructor & Destructor Documentation

27.196.4.1 `gdcm::Parser::Parser ( )` `[inline]`

27.196.4.2 `gdcm::Parser::~~Parser ( )` `[inline]`

### 27.196.5 Member Function Documentation

27.196.5.1 `char* gdcm::Parser::GetBuffer ( int len )` `[protected]`

27.196.5.2 `unsigned long gdcm::Parser::GetCurrentByteIndex ( )` `const`

27.196.5.3 `ErrorType gdcm::Parser::GetErrorCode ( )` `const`

27.196.5.4 `static const char* gdcm::Parser::GetErrorString ( ErrorType const & err )` `[static]`

27.196.5.5 `void* gdcm::Parser::GetUserData ( )` `const`

27.196.5.6 `bool gdcM::Parser::Parse ( const char * s, int len, bool isFinal )`

27.196.5.7 `bool gdcM::Parser::ParseBuffer ( int len, bool isFinal )` [protected]

27.196.5.8 `ErrorType gdcM::Parser::Process ( )` [protected]

27.196.5.9 `void gdcM::Parser::SetElementHandler ( StartElementHandler start, EndElementHandler end )`

27.196.5.10 `void gdcM::Parser::SetUserData ( void * userData )`

The documentation for this class was generated from the following file:

- [gdcMParser.h](#)

## 27.197 gdcM::Patient Class Reference

See PS 3.3 - 2007 DICOM MODEL OF THE REAL-WORLD, p 54.

```
#include <gdcMPatient.h>
```

### Public Member Functions

- [Patient \(\)](#)

#### 27.197.1 Detailed Description

See PS 3.3 - 2007 DICOM MODEL OF THE REAL-WORLD, p 54.

#### 27.197.2 Constructor & Destructor Documentation

27.197.2.1 `gdcM::Patient::Patient ( )` [inline]

The documentation for this class was generated from the following file:

- [gdcMPatient.h](#)


## 27.198 gdcM::network::PDataTFPDU Class Reference

[PDataTFPDU Table](#) 9-22 P-DATA-TF PDU FIELDS.


```
#include <gdcMPDataTFPDU.h>
```



Inheritance diagram for gdcm::network::PDataTFPDU:

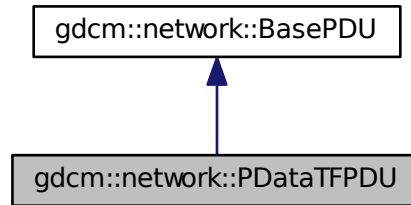

Collaboration diagram for gdcm::network::PDataTFPDU:

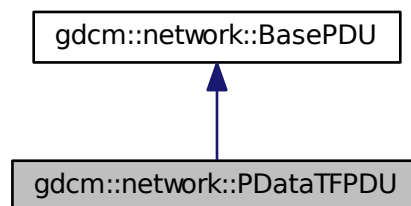

## Public Types

- typedef std::vector< [PresentationDataValue](#) >::size\_type [SizeType](#)

## Public Member Functions

- [PDataTFPDU](#) ()
- void [AddPresentationDataValue](#) ([PresentationDataValue](#) const &pdv)
- [SizeType](#) [GetNumberOfPresentationDataValues](#) () const
- [PresentationDataValue](#) const & [GetPresentationDataValue](#) ([SizeType](#) i) const
- bool [IsLastFragment](#) () const
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

## Protected Member Functions

- `std::istream & ReadInto (std::istream &is, std::ostream &os)`

### 27.198.1 Detailed Description

[PDataTFPDU Table](#) 9-22 P-DATA-TF PDU FIELDS.

### 27.198.2 Member Typedef Documentation

27.198.2.1 `typedef std::vector<PresentationDataValue>::size_type gdcmm::network::PDataTFPDU::SizeType`

### 27.198.3 Constructor & Destructor Documentation

27.198.3.1 `gdcmm::network::PDataTFPDU::PDataTFPDU ( )`

### 27.198.4 Member Function Documentation

27.198.4.1 `void gdcmm::network::PDataTFPDU::AddPresentationDataValue ( PresentationDataValue const & pdv )`  
[inline]

27.198.4.2 `SizeType gdcmm::network::PDataTFPDU::GetNumberOfPresentationDataValues ( ) const` [inline]

27.198.4.3 `PresentationDataValue const& gdcmm::network::PDataTFPDU::GetPresentationDataValue ( SizeType i ) const`  
[inline]

27.198.4.4 `bool gdcmm::network::PDataTFPDU::IsLastFragment ( ) const` [virtual]

Implements [gdcmm::network::BasePDU](#).

27.198.4.5 `void gdcmm::network::PDataTFPDU::Print ( std::ostream & os ) const` [virtual]

Implements [gdcmm::network::BasePDU](#).

27.198.4.6 `std::istream& gdcmm::network::PDataTFPDU::Read ( std::istream & is )` [virtual]

Implements [gdcmm::network::BasePDU](#).

27.198.4.7 `std::istream& gdcmm::network::PDataTFPDU::ReadInto ( std::istream & is, std::ostream & os )` [protected]

27.198.4.8 `size_t gdcmm::network::PDataTFPDU::Size ( ) const` [virtual]

Implements [gdcmm::network::BasePDU](#).

27.198.4.9 `const std::ostream& gdcmm::network::PDataTFPDU::Write ( std::ostream & os ) const` [virtual]

Implements [gdcmm::network::BasePDU](#).

The documentation for this class was generated from the following file:


- [gdcmPDataTFPDU.h](#)

## 27.199 gdcm::PDBElement Class Reference

Class to represent a PDB [Element](#).

```
#include <gdcmPDBElement.h>
```

Collaboration diagram for gdcm::PDBElement:

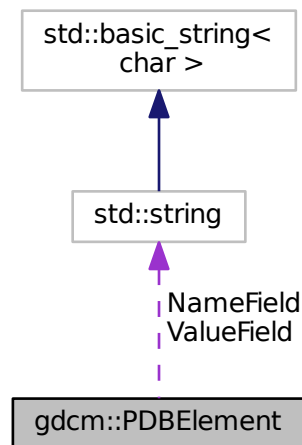

### Public Member Functions

- [PDBElement](#) ()
- const char \* [GetName](#) () const  
*Set/Get Name.*
- const char \* [GetValue](#) () const  
*Set/Get Value.*
- bool [operator==](#) (const [PDBElement](#) &de) const
- void [SetName](#) (const char \*name)
- void [SetValue](#) (const char \*value)

### Protected Attributes

- std::string [NameField](#)
- std::string [ValueField](#)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [PDBElement](#) &val)

### 27.199.1 Detailed Description

Class to represent a PDB [Element](#).

See also

[PDBHeader](#)

### 27.199.2 Constructor & Destructor Documentation

27.199.2.1 `gdcm::PDBElement::PDBElement ( )` `[inline]`

### 27.199.3 Member Function Documentation

27.199.3.1 `const char* gdcm::PDBElement::GetName ( ) const` `[inline]`

Set/Get Name.

27.199.3.2 `const char* gdcm::PDBElement::GetValue ( ) const` `[inline]`

Set/Get [Value](#).

27.199.3.3 `bool gdcm::PDBElement::operator== ( const PDBElement & de ) const` `[inline]`

References NameField, and ValueField.

27.199.3.4 `void gdcm::PDBElement::SetName ( const char * name )` `[inline]`

27.199.3.5 `void gdcm::PDBElement::SetValue ( const char * value )` `[inline]`

### 27.199.4 Friends And Related Function Documentation

27.199.4.1 `std::ostream& operator<< ( std::ostream & os, const PDBElement & val )` `[friend]`

### 27.199.5 Member Data Documentation

27.199.5.1 `std::string gdcm::PDBElement::NameField` `[protected]`

Referenced by `gdcm::operator<<()`, and `operator==()`.

27.199.5.2 `std::string gdcm::PDBElement::ValueField` `[protected]`

Referenced by `gdcm::operator<<()`, and `operator==()`.

The documentation for this class was generated from the following file:

- [gdcmPDBElement.h](#)

## 27.200 gdcm::PDBHeader Class Reference

Class for [PDBHeader](#).

```
#include <gdcmPDBHeader.h>
```

### Public Member Functions

- [PDBHeader](#) ()
- [~PDBHeader](#) ()
- bool [FindPDBElementByName](#) (const char \*name)  
*Return true if the PDB element matching name is found or not.*
- const [PDBElement](#) & [GetPDBElementByName](#) (const char \*name)
- bool [LoadFromDataElement](#) ([DataElement](#) const &de)  
*Load the PDB Header from a [DataElement](#) of a [DataSet](#).*
- void [Print](#) (std::ostream &os) const  
*Print.*

### Static Public Member Functions

- static const [PrivateTag](#) & [GetPDBInfoTag](#) ()  
*Return the Private [Tag](#) where the PDB header is stored within a DICOM [DataSet](#).*

### Protected Member Functions

- const [PDBElement](#) & [GetPDBEEnd](#) () const

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [PDBHeader](#) &d)

### 27.200.1 Detailed Description

Class for [PDBHeader](#).

GEMS MR [Image](#) have an [Attribute](#) (0025,1b,GEMS\_SERS\_01) which store the Acquisition parameter of the MR [Image](#). It is compressed and can therefore not be used as is. This class de-encapsulated the Protocol Data Block and allow users to query element by name.

#### Warning

Everything you do with this code is at your own risk, since decoding process was not written from specification documents.  
: the API of this class might change.

#### See also

[CSAHeader](#)

## 27.200.2 Constructor & Destructor Documentation

27.200.2.1 `gdcm::PDBHeader::PDBHeader ( )` `[inline]`

27.200.2.2 `gdcm::PDBHeader::~~PDBHeader ( )` `[inline]`

## 27.200.3 Member Function Documentation

27.200.3.1 `bool gdcm::PDBHeader::FindPDBElementByName ( const char * name )`

Return true if the PDB element matching name is found or not.

27.200.3.2 `const PDBElement& gdcm::PDBHeader::GetPDBEEnd ( ) const` `[protected]`

27.200.3.3 `const PDBElement& gdcm::PDBHeader::GetPDBElementByName ( const char * name )`

Lookup in the PDB header if a PDB element match the name 'name':

### Warning

Case Sensitive

27.200.3.4 `static const PrivateTag& gdcm::PDBHeader::GetPDBInfoTag ( )` `[static]`

Return the Private [Tag](#) where the PDB header is stored within a DICOM [DataSet](#).

27.200.3.5 `bool gdcm::PDBHeader::LoadFromDataElement ( DataElement const & de )`

Load the PDB Header from a [DataElement](#) of a [DataSet](#).

27.200.3.6 `void gdcm::PDBHeader::Print ( std::ostream & os ) const`

Print.

Referenced by `gdcm::operator<<()`.

## 27.200.4 Friends And Related Function Documentation

27.200.4.1 `std::ostream& operator<< ( std::ostream & _os, const PDBHeader & d )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmPDBHeader.h](#)

## 27.201 gdcm::PDFCodec Class Reference


[PDFCodec](#) class.

```
#include <gdcmPDFCodec.h>
```

Inheritance diagram for gdcm::PDFCodec:

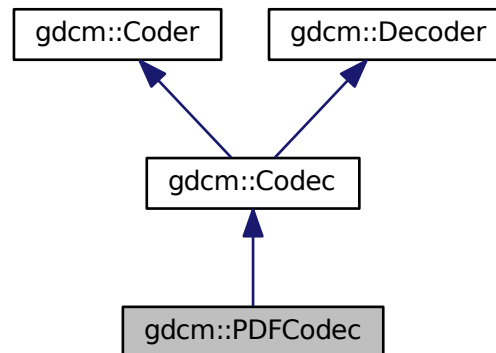

Collaboration diagram for gdcm::PDFCodec:

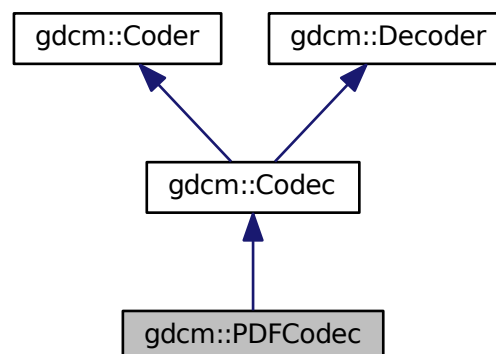

## Public Member Functions

- [PDFCodec](#) ()
- [~PDFCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &) const  
*Return whether this coder support this transfer syntax (can code it)*
- bool [CanDecode](#) ([TransferSyntax](#) const &) const

*Return whether this decoder support this transfer syntax (can decode it)*

- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)

*Decode.*

## Additional Inherited Members

### 27.201.1 Detailed Description

[PDFCodec](#) class.

### 27.201.2 Constructor & Destructor Documentation

27.201.2.1 [gdcm::PDFCodec::PDFCodec](#) ( )

27.201.2.2 [gdcm::PDFCodec::~~PDFCodec](#) ( )

### 27.201.3 Member Function Documentation

27.201.3.1 bool [gdcm::PDFCodec::CanCode](#) ( [TransferSyntax](#) const & ) const [inline],[virtual]

Return whether this coder support this transfer syntax (can code it)

Implements [gdcm::Coder](#).

27.201.3.2 bool [gdcm::PDFCodec::CanDecode](#) ( [TransferSyntax](#) const & ) const [inline],[virtual]

Return whether this decoder support this transfer syntax (can decode it)

Implements [gdcm::Decoder](#).

27.201.3.3 bool [gdcm::PDFCodec::Decode](#) ( [DataElement](#) const & , [DataElement](#) & ) [virtual]

*Decode.*

Reimplemented from [gdcm::Decoder](#).

The documentation for this class was generated from the following file:

- [gdcmPDFCodec.h](#)

## 27.202 gdcm::network::PDUFactory Class Reference

[PDUFactory](#) basically, given an initial byte, construct the appropriate PDU. This way, the event loop doesn't have to know about all the different PDU types.

```
#include <gdcmPDUFactory.h>
```

### Static Public Member Functions

- static [BasePDU](#) \* [ConstructAbortPDU](#) ()



- static [BasePDU](#) \* [ConstructPDU](#) (uint8\_t itemtype)
- static [BasePDU](#) \* [ConstructReleasePDU](#) ()
- static std::vector< [BasePDU](#) \* > [CreateCEchoPDU](#) (const [ULConnection](#) &inConnection)
- static std::vector< [BasePDU](#) \* > [CreateCFindPDU](#) (const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \*inRootQuery)
- static std::vector< [BasePDU](#) \* > [CreateCMovePDU](#) (const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \*inRootQuery)
- static std::vector< [BasePDU](#) \* > [CreateCStoreRQPDU](#) (const [ULConnection](#) &inConnection, const [File](#) &file)
- static std::vector< [BasePDU](#) \* > [CreateCStoreRSPPDU](#) (const [DataSet](#) \*inDataSet, const [BasePDU](#) \*inPC)
- static [EEventID](#) [DetermineEventByPDU](#) (const [BasePDU](#) \*inPDU)
- static std::vector< [PresentationDataValue](#) > [GetPDVs](#) (const std::vector< [BasePDU](#) \* > &inDataPDUs)

### 27.202.1 Detailed Description

[PDUFactory](#) basically, given an initial byte, construct the appropriate PDU. This way, the event loop doesn't have to know about all the different PDU types.

### 27.202.2 Member Function Documentation

- 27.202.2.1 static [BasePDU](#)\* [gdcmm::network::PDUFactory::ConstructAbortPDU](#) ( ) [static]
- 27.202.2.2 static [BasePDU](#)\* [gdcmm::network::PDUFactory::ConstructPDU](#) ( uint8\_t itemtype ) [static]
- 27.202.2.3 static [BasePDU](#)\* [gdcmm::network::PDUFactory::ConstructReleasePDU](#) ( ) [static]
- 27.202.2.4 static std::vector<[BasePDU](#)\*> [gdcmm::network::PDUFactory::CreateCEchoPDU](#) ( const [ULConnection](#) &inConnection ) [static]
- 27.202.2.5 static std::vector<[BasePDU](#)\*> [gdcmm::network::PDUFactory::CreateCFindPDU](#) ( const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \* inRootQuery ) [static]
- 27.202.2.6 static std::vector<[BasePDU](#)\*> [gdcmm::network::PDUFactory::CreateCMovePDU](#) ( const [ULConnection](#) &inConnection, const [BaseRootQuery](#) \* inRootQuery ) [static]
- 27.202.2.7 static std::vector<[BasePDU](#)\*> [gdcmm::network::PDUFactory::CreateCStoreRQPDU](#) ( const [ULConnection](#) &inConnection, const [File](#) & file ) [static]
- 27.202.2.8 static std::vector<[BasePDU](#)\*> [gdcmm::network::PDUFactory::CreateCStoreRSPPDU](#) ( const [DataSet](#) \* inDataSet, const [BasePDU](#) \* inPC ) [static]
- 27.202.2.9 static [EEventID](#) [gdcmm::network::PDUFactory::DetermineEventByPDU](#) ( const [BasePDU](#) \* inPDU ) [static]
- 27.202.2.10 static std::vector<[PresentationDataValue](#)> [gdcmm::network::PDUFactory::GetPDVs](#) ( const std::vector< [BasePDU](#) \* > &inDataPDUs ) [static]

The documentation for this class was generated from the following file:

- [gdcmmPDUFactory.h](#)

## 27.203 gdcm::PersonName Class Reference

[PersonName](#) class.

```
#include <gdcmPersonName.h>
```

### Public Member Functions

- unsigned int [GetMaxLength](#) () const
- unsigned int [GetNumberOfComponents](#) () const
- void [Print](#) (std::ostream &os) const
- void [SetBlob](#) (const std::vector< char > &v)
- void [SetComponents](#) (const char \*comp1="", const char \*comp2="", const char \*comp3="", const char \*comp4="", const char \*comp5="")
- void [SetComponents](#) (const char \*components[])

### Public Attributes

- char [Component](#) [[MaxNumberOfComponents](#)][[MaxLength](#)+1]

### Static Public Attributes

- static const unsigned int [MaxLength](#) = 64
- static const unsigned int [MaxNumberOfComponents](#) = 5
- static const char [Padding](#) = ' '
- static const char [Separator](#) = '^'

### 27.203.1 Detailed Description

[PersonName](#) class.

### 27.203.2 Member Function Documentation

27.203.2.1 unsigned int gdcm::PersonName::GetMaxLength ( ) const [\[inline\]](#)

27.203.2.2 unsigned int gdcm::PersonName::GetNumberOfComponents ( ) const [\[inline\]](#)

27.203.2.3 void gdcm::PersonName::Print ( std::ostream & os ) const [\[inline\]](#)

27.203.2.4 void gdcm::PersonName::SetBlob ( const std::vector< char > & v ) [\[inline\]](#)

27.203.2.5 void gdcm::PersonName::SetComponents ( const char \* *comp1* = " ", const char \* *comp2* = " ", const char \* *comp3* = " ", const char \* *comp4* = " ", const char \* *comp5* = " " ) [\[inline\]](#)

27.203.2.6 void gdcm::PersonName::SetComponents ( const char \* *components*[] ) [\[inline\]](#)

### 27.203.3 Member Data Documentation

27.203.3.1 `char gdcm::PersonName::Component[MaxNumberOfComponents][MaxLength+1]`

27.203.3.2 `const unsigned int gdcm::PersonName::MaxLength = 64` `[static]`

27.203.3.3 `const unsigned int gdcm::PersonName::MaxNumberOfComponents = 5` `[static]`

27.203.3.4 `const char gdcm::PersonName::Padding = ''` `[static]`

27.203.3.5 `const char gdcm::PersonName::Separator = '^'` `[static]`

The documentation for this class was generated from the following file:


- [gdcmPersonName.h](#)

## 27.204 gdcm::PGXCodec Class Reference


Class to do PGX See PGX as used in JPEG 2000 implementation and reference images.

```
#include <gdcmPGXCodec.h>
```

Inheritance diagram for gdcm::PGXCodec:

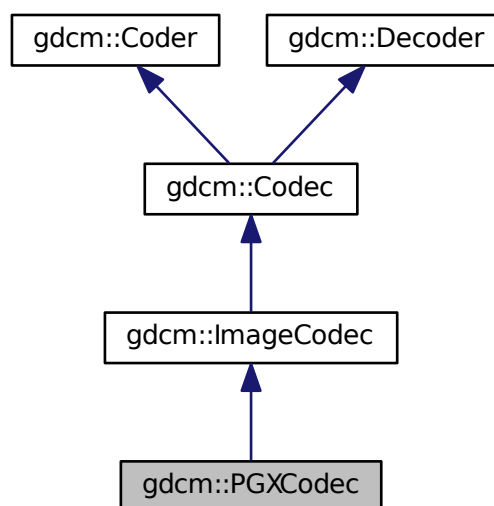

Collaboration diagram for `gdcm::PGXCodec`:

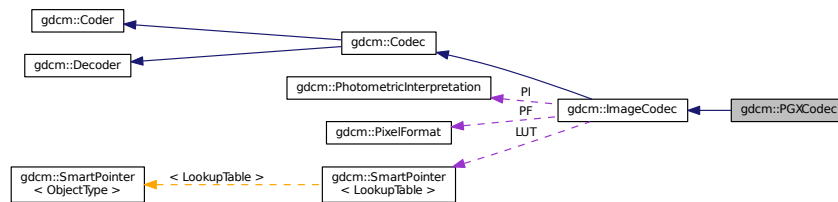

## Public Member Functions

- `PGXCodec ()`
- `~PGXCodec ()`
- `bool CanCode (TransferSyntax const &ts) const`  
*Return whether this coder support this transfer syntax (can code it)*
- `bool CanDecode (TransferSyntax const &ts) const`  
*Return whether this decoder support this transfer syntax (can decode it)*
- `virtual ImageCodec * Clone () const`
- `bool GetHeaderInfo (std::istream &is, TransferSyntax &ts)`
- `bool Read (const char *filename, DataElement &out) const`
- `bool Write (const char *filename, const DataElement &out) const`

## Additional Inherited Members

### 27.204.1 Detailed Description

Class to do PGX See PGX as used in JPEG 2000 implementation and reference images.

### 27.204.2 Constructor & Destructor Documentation

27.204.2.1 `gdcm::PGXCodec::PGXCodec ( )`

27.204.2.2 `gdcm::PGXCodec::~~PGXCodec ( )`

### 27.204.3 Member Function Documentation

27.204.3.1 `bool gdcm::PGXCodec::CanCode ( TransferSyntax const & ) const` `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Reimplemented from `gdcm::ImageCodec`.

27.204.3.2 `bool gdcm::PGXCodec::CanDecode ( TransferSyntax const & ) const` `[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from `gdcm::ImageCodec`.

27.204.3.3 virtual ImageCodec\* gdcm::PGXCodec::Clone ( ) const [virtual]

Implements [gdcm::ImageCodec](#).

27.204.3.4 bool gdcm::PGXCodec::GetHeaderInfo ( std::istream & *is*, TransferSyntax & *ts* ) [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.204.3.5 bool gdcm::PGXCodec::Read ( const char \* *filename*, DataElement & *out* ) const

27.204.3.6 bool gdcm::PGXCodec::Write ( const char \* *filename*, const DataElement & *out* ) const

The documentation for this class was generated from the following file:

- [gdcmPGXCodec.h](#)

## 27.205 gdcm::PhotometricInterpretation Class Reference

Class to represent an [PhotometricInterpretation](#).

```
#include <gdcmPhotometricInterpretation.h>
```

### Public Types

- enum [PType](#) {  
[UNKNOWN](#) = 0,  
[MONOCHROME1](#),  
[MONOCHROME2](#),  
[PALETTE\\_COLOR](#),  
[RGB](#),  
[HSV](#),  
[ARGB](#),  
[CMYK](#),  
[YBR\\_FULL](#),  
[YBR\\_FULL\\_422](#),  
[YBR\\_PARTIAL\\_422](#),  
[YBR\\_PARTIAL\\_420](#),  
[YBR\\_ICT](#),  
[YBR\\_RCT](#),  
[PI\\_END](#) }

### Public Member Functions

- [PhotometricInterpretation](#) ([PType](#) *pi*=[UNKNOWN](#))
- unsigned short [GetSamplesPerPixel](#) () const  
*return the value for Sample Per Pixel associated with a particular Photometric Interpretation*
- const char \* [GetString](#) () const
- [PType](#) [GetType](#) () const
- bool [IsLossless](#) () const

- bool [IsLossy](#) () const
- bool [IsSameColorSpace](#) ([PhotometricInterpretation](#) const &pi) const
- [operator PType](#) () const

### Static Public Member Functions

- static const char \* [GetPIString](#) ([PType](#) pi)
- static [PType](#) [GetPType](#) (const char \*pi)
- static bool [IsRetired](#) ([PType](#) pi)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [PhotometricInterpretation](#) &pi)

## 27.205.1 Detailed Description

Class to represent an [PhotometricInterpretation](#).

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [HelloVizWorld.cxx](#), and [iU22tomultisc.cxx](#).

## 27.205.2 Member Enumeration Documentation

### 27.205.2.1 enum gdcm::PhotometricInterpretation::PType

Enumerator

***UNKNOWN***  
***MONOCHROME1***  
***MONOCHROME2***  
***PALETTE\_COLOR***  
***RGB***  
***HSV***  
***ARGB***  
***CMYK***  
***YBR\_FULL***  
***YBR\_FULL\_422***  
***YBR\_PARTIAL\_422***  
***YBR\_PARTIAL\_420***  
***YBR\_ICT***  
***YBR\_RCT***  
***PI\_END***

### 27.205.3 Constructor & Destructor Documentation

27.205.3.1 `gdcm::PhotometricInterpretation::PhotometricInterpretation ( PType pi = UNKNOWN )` `[inline]`

### 27.205.4 Member Function Documentation

27.205.4.1 `static const char* gdcm::PhotometricInterpretation::GetPIString ( PType pi )` `[static]`

Referenced by `gdcm::operator<<()`.

27.205.4.2 `static PType gdcm::PhotometricInterpretation::GetPType ( const char * pi )` `[static]`

27.205.4.3 `unsigned short gdcm::PhotometricInterpretation::GetSamplesPerPixel ( )` `const`

return the value for Sample Per Pixel associated with a particular Photometric Interpretation

27.205.4.4 `const char* gdcm::PhotometricInterpretation::GetString ( )` `const`

27.205.4.5 `PType gdcm::PhotometricInterpretation::GetType ( )` `const` `[inline]`

27.205.4.6 `bool gdcm::PhotometricInterpretation::IsLossless ( )` `const`

27.205.4.7 `bool gdcm::PhotometricInterpretation::IsLossy ( )` `const`

27.205.4.8 `static bool gdcm::PhotometricInterpretation::IsRetired ( PType pi )` `[static]`

27.205.4.9 `bool gdcm::PhotometricInterpretation::IsSameColorSpace ( PhotometricInterpretation const & pi )` `const`

27.205.4.10 `gdcm::PhotometricInterpretation::operator PType ( )` `const` `[inline]`

### 27.205.5 Friends And Related Function Documentation

27.205.5.1 `std::ostream& operator<< ( std::ostream & os, const PhotometricInterpretation & pi )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmPhotometricInterpretation.h](#)

## 27.206 gdcm::PixelFormat Class Reference

[PixelFormat](#).

```
#include <gdcmPixelFormat.h>
```

## Public Types

- enum [ScalarType](#) {  
[UINT8](#),  
[INT8](#),  
[UINT12](#),  
[INT12](#),  
[UINT16](#),  
[INT16](#),  
[UINT32](#),  
[INT32](#),  
[UINT64](#),  
[INT64](#),  
[FLOAT16](#),  
[FLOAT32](#),  
[FLOAT64](#),  
[SINGLEBIT](#),  
[UNKNOWN](#) }

## Public Member Functions

- [PixelFormat](#) (unsigned short samplesperpixel=1, unsigned short bitsallocated=8, unsigned short bitsstored=8, unsigned short highbit=7, unsigned short pixelrepresentation=0)
- [PixelFormat](#) ([ScalarType](#) st)
- unsigned short [GetBitsAllocated](#) () const  
*BitsAllocated see [Tag](#) (0028,0100) US Bits Allocated.*
- unsigned short [GetBitsStored](#) () const  
*BitsStored see [Tag](#) (0028,0101) US Bits Stored.*
- unsigned short [GetHighBit](#) () const  
*HighBit see [Tag](#) (0028,0102) US High Bit.*
- int64\_t [GetMax](#) () const  
*return the max possible of the pixel*
- int64\_t [GetMin](#) () const  
*return the min possible of the pixel*
- unsigned short [GetPixelRepresentation](#) () const  
*PixelRepresentation: 0 or 1, see [Tag](#) (0028,0103) US Pixel Representation.*
- uint8\_t [GetPixelSize](#) () const
- unsigned short [GetSamplesPerPixel](#) () const
- [ScalarType](#) [GetScalarType](#) () const  
*ScalarType does not take into account the sample per pixel.*
- const char \* [GetScalarTypeAsString](#) () const
- bool [IsCompatible](#) (const [TransferSyntax](#) &ts) const
- bool [IsValid](#) () const  
*return IsValid*
- operator [ScalarType](#) () const
- bool operator!= ([ScalarType](#) st) const
- bool operator!= (const [PixelFormat](#) &pf) const
- bool operator== ([ScalarType](#) st) const
- bool operator== (const [PixelFormat](#) &pf) const
- void [Print](#) (std::ostream &os) const



*Print.*

- void [SetBitsAllocated](#) (unsigned short ba)
- void [SetBitsStored](#) (unsigned short bs)
- void [SetHighBit](#) (unsigned short hb)
- void [SetPixelRepresentation](#) (unsigned short pr)
- void [SetSamplesPerPixel](#) (unsigned short spp)
- void [SetScalarType](#) ([ScalarType](#) st)

## Protected Member Functions

- bool [Validate](#) ()

*When image with 24/24/23 was read, need to validate.*

## Friends

- class [Bitmap](#)
- std::ostream & [operator<<](#) (std::ostream &\_os, const [PixelFormat](#) &pf)

## 27.206.1 Detailed Description

[PixelFormat](#).

### Note

By default the Pixel [Type](#) will be instantiated with the following parameters:

- SamplesPerPixel : 1
- BitsAllocated : 8
- BitsStored : 8
- HighBit : 7
- PixelRepresentation : 0

Fundamentally [PixelFormat](#) is very close to what DICOM allows. It will be very hard to extend this class for the upcoming DICOM standard where Floating 32 and 64bits will be allowed.

It is also very hard for this class to fully support 64bits integer type (see GetMin / GetMax signature restricted to 64bits signed).

### Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [csa2img.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GetJPEGSample←Precision.cxx](#), [iU22tomultisc.cxx](#), and [threadgdcm.cxx](#).

## 27.206.2 Member Enumeration Documentation

### 27.206.2.1 enum gdcm::PixelFormat::ScalarType

#### Enumerator

**UINT8**

**INT8**

**UINT12**  
**INT12**  
**UINT16**  
**INT16**  
**UINT32**  
**INT32**  
**UINT64**  
**INT64**  
**FLOAT16**  
**FLOAT32**  
**FLOAT64**  
**SINGLEBIT**  
**UNKNOWN**

### 27.206.3 Constructor & Destructor Documentation

27.206.3.1 `gdcm::PixelFormat::PixelFormat ( unsigned short samplesperpixel = 1, unsigned short bitsallocated = 8, unsigned short bitsstored = 8, unsigned short highbit = 7, unsigned short pixelrepresentation = 0 )` `[inline]`, `[explicit]`

27.206.3.2 `gdcm::PixelFormat::PixelFormat ( ScalarType st )`

### 27.206.4 Member Function Documentation

27.206.4.1 `unsigned short gdcm::PixelFormat::GetBitsAllocated ( ) const` `[inline]`

BitsAllocated see [Tag](#) (0028,0100) US Bits Allocated.

Examples:

[GetJPEGSamplePrecision.cxx](#).

27.206.4.2 `unsigned short gdcm::PixelFormat::GetBitsStored ( ) const` `[inline]`

BitsStored see [Tag](#) (0028,0101) US Bits Stored.

Examples:

[GetJPEGSamplePrecision.cxx](#).

27.206.4.3 `unsigned short gdcm::PixelFormat::GetHighBit ( ) const` `[inline]`

HighBit see [Tag](#) (0028,0102) US High Bit.

27.206.4.4 `int64_t gdcm::PixelFormat::GetMax ( ) const`

return the max possible of the pixel

27.206.4.5 `int64_t gdcm::PixelFormat::GetMin ( ) const`

return the min possible of the pixel

27.206.4.6 `unsigned short gdcm::PixelFormat::GetPixelRepresentation ( ) const` `[inline]`

PixelRepresentation: 0 or 1, see [Tag](#) (0028,0103) US Pixel Representation.

27.206.4.7 `uint8_t gdcm::PixelFormat::GetPixelSize ( ) const`

return the size of the pixel This is the number of words it would take to store one pixel

#### Warning

the return value takes into account the SamplesPerPixel  
in the rare case when BitsAllocated == 12, the function assume word padding and value returned will be identical  
as if BitsAllocated == 16

#### Examples:

[threadgdcm.cxx](#).

27.206.4.8 `unsigned short gdcm::PixelFormat::GetSamplesPerPixel ( ) const`

Samples Per Pixel see (0028,0002) US Samples Per Pixel DICOM - only allows 1, 3 and 4 as valid value. Other value are undefined behavior.

#### Examples:

[threadgdcm.cxx](#).

27.206.4.9 `ScalarType gdcm::PixelFormat::GetScalarType ( ) const`

ScalarType does not take into account the sample per pixel.

27.206.4.10 `const char* gdcm::PixelFormat::GetScalarTypeAsString ( ) const`

27.206.4.11 `bool gdcm::PixelFormat::IsCompatible ( const TransferSyntax & ts ) const`

27.206.4.12 `bool gdcm::PixelFormat::IsValid ( ) const`

return IsValid

27.206.4.13 `gdcm::PixelFormat::operator ScalarType ( ) const` `[inline]`

27.206.4.14 `bool gdcm::PixelFormat::operator!= ( ScalarType st ) const` `[inline]`

27.206.4.15 `bool gdcm::PixelFormat::operator!= ( const PixelFormat & pf ) const` `[inline]`

27.206.4.16 `bool gdcm::PixelFormat::operator==( ScalarType st ) const` `[inline]`

27.206.4.17 `bool gdcm::PixelFormat::operator==( const PixelFormat & pf ) const` `[inline]`

27.206.4.18 `void gdcm::PixelFormat::Print ( std::ostream & os ) const`

Print.

Referenced by `gdcm::operator<<()`.

27.206.4.19 `void gdcm::PixelFormat::SetBitsAllocated ( unsigned short ba )` `[inline]`

27.206.4.20 `void gdcm::PixelFormat::SetBitsStored ( unsigned short bs )` `[inline]`

27.206.4.21 `void gdcm::PixelFormat::SetHighBit ( unsigned short hb )` `[inline]`

27.206.4.22 `void gdcm::PixelFormat::SetPixelRepresentation ( unsigned short pr )` `[inline]`

27.206.4.23 `void gdcm::PixelFormat::SetSamplesPerPixel ( unsigned short spp )` `[inline]`

Examples:

[CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), and [GenFakelImage.cxx](#).

References `gdcmAssertMacro`.

27.206.4.24 `void gdcm::PixelFormat::SetScalarType ( ScalarType st )`

Set [PixelFormat](#) based only on the `ScalarType`

Warning

: You need to call `SetScalarType` *before* `SetSamplesPerPixel`

27.206.4.25 `bool gdcm::PixelFormat::Validate ( )` `[protected]`

When image with 24/24/23 was read, need to validate.

Referenced by `gdcm::Bitmap::SetPixelFormat()`.

## 27.206.5 Friends And Related Function Documentation

27.206.5.1 `friend class Bitmap` `[friend]`


27.206.5.2 `std::ostream& operator<< ( std::ostream & _os, const PixelFormat & pf )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmPixelFormat.h](#)

**Pixmap** class A bitmap based image. Used as parent for both IconImage and the main Pixel Data **Image** It does not contains any World Space information (IPP, IOP)

Inheritance diagram for `gdcm::Pixmap`:

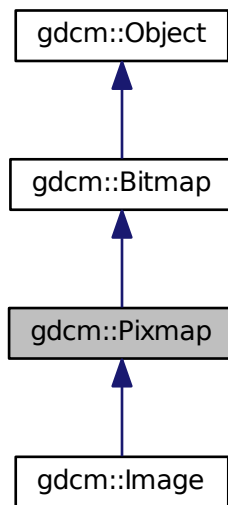[illegible]

- Pixmap ()

- [~Pixmap](#) ()
- bool [AreOverlaysInPixelData](#) () const  
*returns if Overlays are stored in the unused bit of the pixel data:*
- [Curve](#) & [GetCurve](#) (size\_t i=0)  
*Curve: group 50xx.*
- const [Curve](#) & [GetCurve](#) (size\_t i=0) const
- const [IconImage](#) & [GetIconImage](#) () const  
*Set/Get Icon Image.*
- [IconImage](#) & [GetIconImage](#) ()
- size\_t [GetNumberOfCurves](#) () const
- size\_t [GetNumberOfOverlays](#) () const
- [Overlay](#) & [GetOverlay](#) (size\_t i=0)  
*Overlay: group 60xx.*
- const [Overlay](#) & [GetOverlay](#) (size\_t i=0) const
- void [Print](#) (std::ostream &) const
- void [RemoveOverlay](#) (size\_t i)
- void [SetIconImage](#) ([IconImage](#) const &ii)
- void [SetNumberOfCurves](#) (size\_t n)
- void [SetNumberOfOverlays](#) (size\_t n)

## Protected Attributes

- std::vector< [Curve](#) > [Curves](#)
- [SmartPointer](#)< [IconImage](#) > [Icon](#)
- std::vector< [Overlay](#) > [Overlays](#)

## Additional Inherited Members

### 27.207.1 Detailed Description

[Pixmap](#) class A bitmap based image. Used as parent for both [IconImage](#) and the main Pixel Data [Image](#) It does not contains any World Space information (IPP, IOP)

See also

[PixmapReader](#)

### 27.207.2 Constructor & Destructor Documentation

27.207.2.1 [gdcm::Pixmap::Pixmap](#) ( )

27.207.2.2 [gdcm::Pixmap::~~Pixmap](#) ( )

### 27.207.3 Member Function Documentation

27.207.3.1 bool [gdcm::Pixmap::AreOverlaysInPixelData](#) ( ) const [virtual]

returns if Overlays are stored in the unused bit of the pixel data:

Reimplemented from [gdcm::Bitmap](#).

27.207.3.2 **Curve&** gdcm::Pixmap::GetCurve ( size\_t *i* = 0 ) [inline]

[Curve](#): group 50xx.

27.207.3.3 **const Curve&** gdcm::Pixmap::GetCurve ( size\_t *i* = 0 ) **const** [inline]

27.207.3.4 **const IconImage&** gdcm::Pixmap::GetIconImage ( ) **const** [inline]

Set/Get Icon [Image](#).

27.207.3.5 **IconImage&** gdcm::Pixmap::GetIconImage ( ) [inline]

27.207.3.6 **size\_t** gdcm::Pixmap::GetNumberOfCurves ( ) **const** [inline]

27.207.3.7 **size\_t** gdcm::Pixmap::GetNumberOfOverlays ( ) **const** [inline]

27.207.3.8 **Overlay&** gdcm::Pixmap::GetOverlay ( size\_t *i* = 0 ) [inline]

[Overlay](#): group 60xx.

27.207.3.9 **const Overlay&** gdcm::Pixmap::GetOverlay ( size\_t *i* = 0 ) **const** [inline]

27.207.3.10 **void** gdcm::Pixmap::Print ( std::ostream & ) **const** [virtual]

Reimplemented from [gdcm::Bitmap](#).

27.207.3.11 **void** gdcm::Pixmap::RemoveOverlay ( size\_t *i* ) [inline]

27.207.3.12 **void** gdcm::Pixmap::SetIconImage ( IconImage **const** & *ii* ) [inline]

27.207.3.13 **void** gdcm::Pixmap::SetNumberOfCurves ( size\_t *n* ) [inline]

27.207.3.14 **void** gdcm::Pixmap::SetNumberOfOverlays ( size\_t *n* ) [inline]

## 27.207.4 Member Data Documentation

27.207.4.1 **std::vector<Curve>** gdcm::Pixmap::Curves [protected]

27.207.4.2 **SmartPointer<IconImage>** gdcm::Pixmap::Icon [protected]

27.207.4.3 **std::vector<Overlay>** gdcm::Pixmap::Overlays [protected]

The documentation for this class was generated from the following file:


- [gdcmPixmap.h](#)

## 27.208 gdcm::PixmapReader Class Reference

[PixmapReader](#).


```
#include <gdcmPixmapReader.h>
```

Inheritance diagram for `gdcm::PixmapReader`:

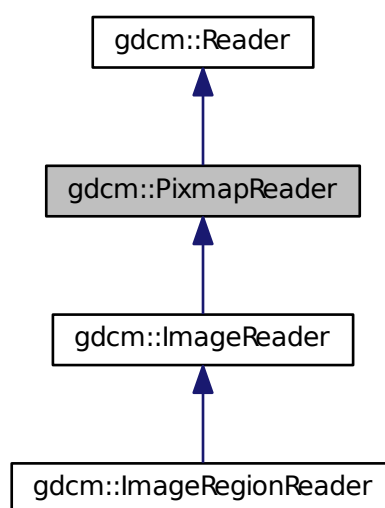



Collaboration diagram for gdcm::PixmapReader:

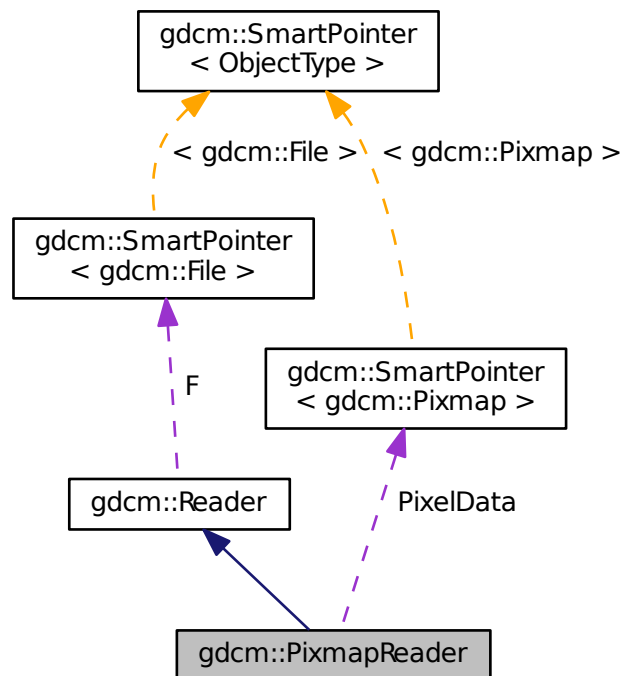

## Public Member Functions

- [PixmapReader](#) ()
- virtual [~PixmapReader](#) ()
- const [Pixmap](#) & [GetPixmap](#) () const  
*Return the read image (need to call [Read\(\)](#) first)*
- [Pixmap](#) & [GetPixmap](#) ()
- virtual bool [Read](#) ()

## Protected Member Functions

- virtual bool [ReadACRNEMAIImage](#) ()
- virtual bool [ReadImage](#) ([MediaStorage](#) const &ms)
- bool [ReadImageInternal](#) ([MediaStorage](#) const &ms, bool handlepixeldata=true)

## Protected Attributes

- [SmartPointer](#) < [Pixmap](#) > [PixelData](#)

### 27.208.1 Detailed Description

[PixmapReader](#).

#### Note

its role is to convert the DICOM [DataSet](#) into a [Pixmap](#) representation By default it is also loading the lookup table and overlay when found as they impact the rendering of the image

See PS 3.3-2008, [Table C.7-11b](#) IMAGE PIXEL MACRO ATTRIBUTES for the list of attribute that belong to what gdcm calls a 'Pixmap'

#### Warning

the API `ReadUpToTag` and `ReadSelectedTag`

#### See also

[Pixmap](#)

### 27.208.2 Constructor & Destructor Documentation

27.208.2.1 `gdcm::PixmapReader::PixmapReader ( )`

27.208.2.2 `virtual gdcm::PixmapReader::~~PixmapReader ( ) [virtual]`

### 27.208.3 Member Function Documentation

27.208.3.1 `const Pixmap& gdcm::PixmapReader::GetPixmap ( ) const`

Return the read image (need to call [Read\(\)](#) first)

27.208.3.2 `Pixmap& gdcm::PixmapReader::GetPixmap ( )`

27.208.3.3 `virtual bool gdcm::PixmapReader::Read ( ) [virtual]`

Read the DICOM image. There are two reason for failure:

1. The input filename is not DICOM
2. The input DICOM file does not contains an [Pixmap](#).

Reimplemented from [gdcm::Reader](#).

Reimplemented in [gdcm::ImageRegionReader](#), and [gdcm::ImageReader](#).

27.208.3.4 `virtual bool gdcm::PixmapReader::ReadACRNEMAIImage ( ) [protected], [virtual]`

Reimplemented in [gdcm::ImageReader](#).

27.208.3.5 `virtual bool gdcm::PixmapReader::ReadImage ( MediaStorage const & ms ) [protected], [virtual]`

Reimplemented in [gdcm::ImageReader](#).

27.208.3.6 `bool gdcm::PixmapReader::ReadImageInternal ( MediaStorage const & ms, bool handlepixeldata = true )`  
`[protected]`

## 27.208.4 Member Data Documentation

27.208.4.1 `SmartPointer<Pixmap> gdcm::PixmapReader::PixelData` `[protected]`

The documentation for this class was generated from the following file:


- [gdcmPixmapReader.h](#)

## 27.209 gdcm::PixmapToPixmapFilter Class Reference


[PixmapToPixmapFilter](#) class Super class for all filter taking an image and producing an output image.

```
#include <gdcmPixmapToPixmapFilter.h>
```

Inheritance diagram for `gdcm::PixmapToPixmapFilter`:

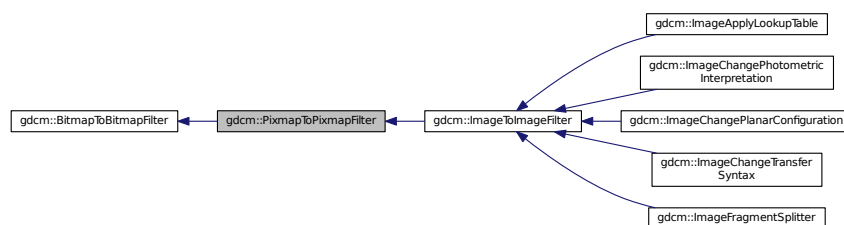

Collaboration diagram for `gdcm::PixmapToPixmapFilter`:

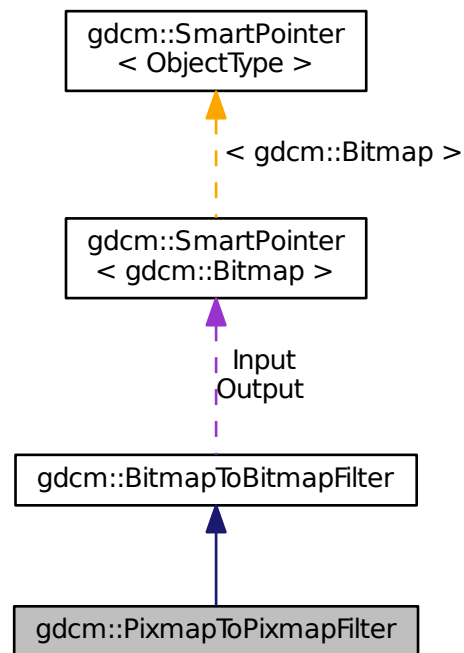

## Public Member Functions

- `PixmapToPixmapFilter ()`
- `~PixmapToPixmapFilter ()`
- `Pixmap & GetInput ()`
- `const Pixmap & GetOutput () const`  
*Get Output image.*
- `const Pixmap & GetOutputAsPixmap () const`

## Additional Inherited Members

### 27.209.1 Detailed Description

`PixmapToPixmapFilter` class Super class for all filter taking an image and producing an output image.

### 27.209.2 Constructor & Destructor Documentation

#### 27.209.2.1 `gdcm::PixmapToPixmapFilter::PixmapToPixmapFilter ( )`

27.209.2.2 `gdcm::PixmapToPixmapFilter::~~PixmapToPixmapFilter ( ) [inline]`

### 27.209.3 Member Function Documentation

27.209.3.1 `Pixmap& gdcm::PixmapToPixmapFilter::GetInput ( )`

27.209.3.2 `const Pixmap& gdcm::PixmapToPixmapFilter::GetOutput ( ) const`

Get Output image.

27.209.3.3 `const Pixmap& gdcm::PixmapToPixmapFilter::GetOutputAsPixmap ( ) const`

The documentation for this class was generated from the following file:


- [gdcmPixmapToPixmapFilter.h](#)

## 27.210 gdcm::PixmapWriter Class Reference


[PixmapWriter](#) This class will takes two inputs:

```
#include <gdcmPixmapWriter.h>
```

Inheritance diagram for `gdcm::PixmapWriter`:

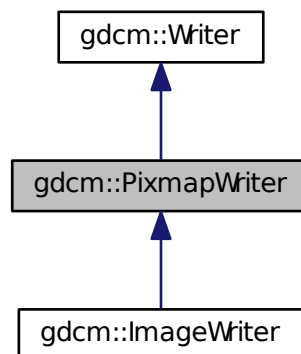

Collaboration diagram for `gdcm::PixmapWriter`:

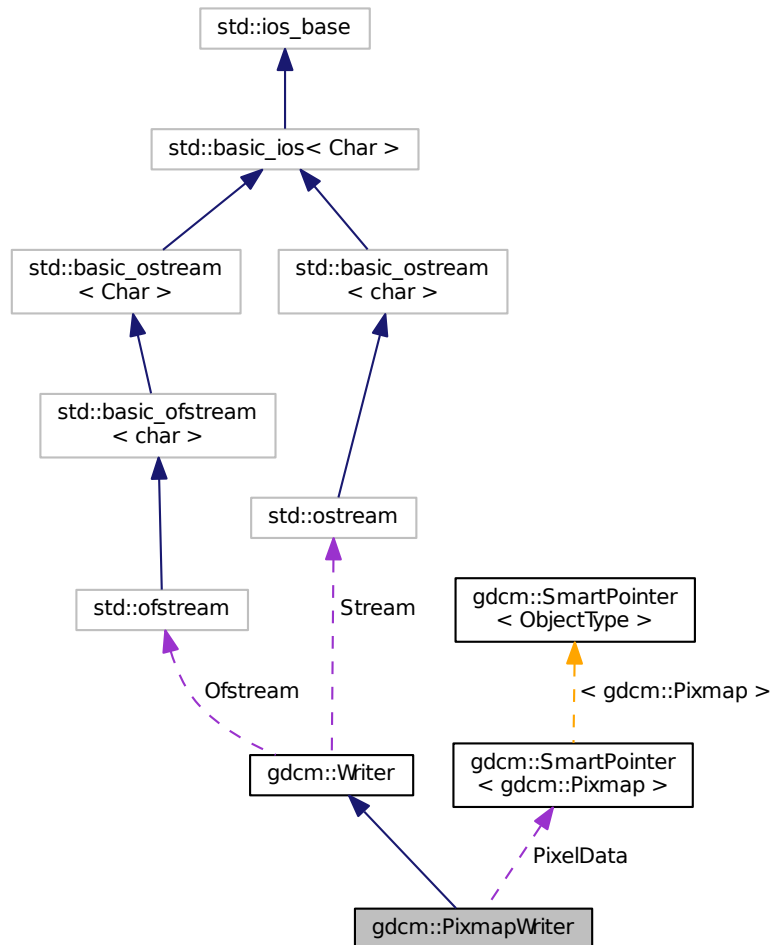

## Public Member Functions

- `PixmapWriter ()`
- `~PixmapWriter ()`
- virtual const `Pixmap & GetImage () const`
- virtual `Pixmap & GetImage ()`
- const `Pixmap & GetPixmap () const`
- `Pixmap & GetPixmap ()`
- virtual void `SetImage (Pixmap const &img)`
- void `SetPixmap (Pixmap const &img)`
- bool `Write ()`

*Write.*

## Protected Member Functions

- void [DolconImage](#) ([DataSet](#) &ds, [Pixmap](#) const &image)
- bool [PrepareWrite](#) ()

## Protected Attributes

- [SmartPointer](#)< [Pixmap](#) > [PixelData](#)

### 27.210.1 Detailed Description

[PixmapWriter](#) This class will takes two inputs:

1. The DICOM [DataSet](#)
2. The [Image](#) input It will override any info from the [Image](#) over the [DataSet](#).

For instance when one read in a lossy compressed image and write out as unencapsulated (ie implicitly lossless) then some attribute are definitely needed to mark this dataset as Lossy (typically 0028,2114)

### 27.210.2 Constructor & Destructor Documentation

27.210.2.1 `gdcm::PixmapWriter::PixmapWriter ( )`

27.210.2.2 `gdcm::PixmapWriter::~~PixmapWriter ( )`

### 27.210.3 Member Function Documentation

27.210.3.1 `void gdcm::PixmapWriter::DolconImage ( DataSet & ds, Pixmap const & image )` `[protected]`

27.210.3.2 `virtual const Pixmap& gdcm::PixmapWriter::GetImage ( ) const` `[inline], [virtual]`

Set/Get [Pixmap](#) to be written It will overwrite anything [Pixmap](#) infos found in [DataSet](#) (see parent class to see how to pass dataset)

Reimplemented in [gdcm::ImageWriter](#).

27.210.3.3 `virtual Pixmap& gdcm::PixmapWriter::GetImage ( )` `[inline], [virtual]`

Reimplemented in [gdcm::ImageWriter](#).

27.210.3.4 `const Pixmap& gdcm::PixmapWriter::GetPixmap ( ) const` `[inline]`

27.210.3.5 `Pixmap& gdcm::PixmapWriter::GetPixmap ( )` `[inline]`

27.210.3.6 `bool gdcm::PixmapWriter::PrepareWrite ( )` `[protected]`

27.210.3.7 virtual void gdcmm::PixmapWriter::SetImage ( Pixmap const & *img* ) [virtual]

Examples:

[CompressImage.cxx](#), [GenFakeImage.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), and [MergeTwoFiles.cxx](#).

27.210.3.8 void gdcmm::PixmapWriter::SetPixmap ( Pixmap const & *img* )

27.210.3.9 bool gdcmm::PixmapWriter::Write ( ) [virtual]

Write.

Reimplemented from [gdcmm::Writer](#).

## 27.210.4 Member Data Documentation

27.210.4.1 SmartPointer<Pixmap> gdcmm::PixmapWriter::PixelData [protected]

The documentation for this class was generated from the following file:

- [gdcmmPixmapWriter.h](#)


## 27.211 gdcmm::PNMCodec Class Reference

Class to do PNM PNM is the Portable anymap file format. The main web page can be found at: <http://netpbm.sourceforge.net/>.↵


```
#include <gdcmmPNMCodec.h>
```



Inheritance diagram for gdcm::PNMCodec:

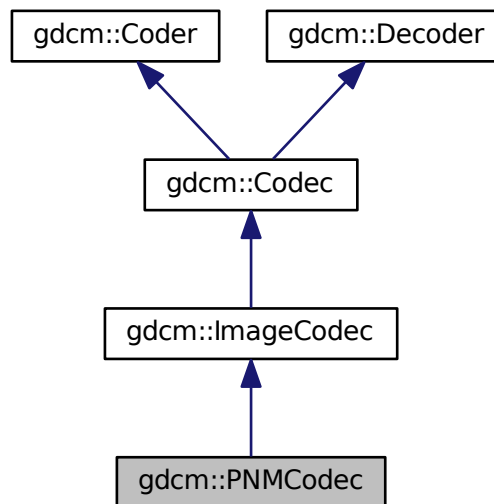

Collaboration diagram for gdcm::PNMCodec:

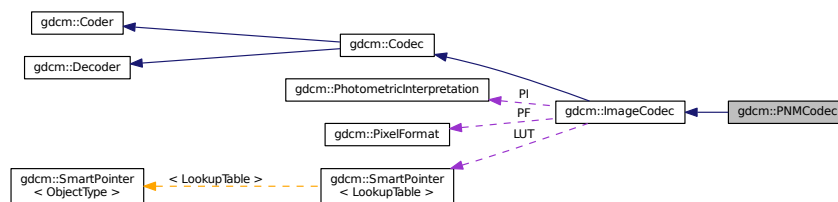

## Public Member Functions

- [PNMCodec](#) ()
- [~PNMCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &ts) const  
Return whether this coder support this transfer syntax (can code it)
- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const  
Return whether this decoder support this transfer syntax (can decode it)
- virtual [ImageCodec](#) \* [Clone](#) () const
- unsigned long [GetBufferLength](#) () const
- bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- bool [Read](#) (const char \*filename, [DataElement](#) &out) const
- void [SetBufferLength](#) (unsigned long l)
- bool [Write](#) (const char \*filename, const [DataElement](#) &out) const

## Additional Inherited Members

### 27.211.1 Detailed Description

Class to do PNM PNM is the Portable anymap file format. The main web page can be found at: <http://netpbm.sourceforge.net/>.

#### Note

Only support P5 & P6 PNM file (binary grayscale and binary rgb)

#### Examples:

[ExtractIconFromFile.cxx](#).

### 27.211.2 Constructor & Destructor Documentation

27.211.2.1 `gdcm::PNMCodec::PNMCodec ( )`

27.211.2.2 `gdcm::PNMCodec::~~PNMCodec ( )`

### 27.211.3 Member Function Documentation

27.211.3.1 `bool gdcm::PNMCodec::CanCode ( TransferSyntax const & ) const` `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.211.3.2 `bool gdcm::PNMCodec::CanDecode ( TransferSyntax const & ) const` `[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.211.3.3 `virtual ImageCodec* gdcm::PNMCodec::Clone ( ) const` `[virtual]`

Implements [gdcm::ImageCodec](#).

27.211.3.4 `unsigned long gdcm::PNMCodec::GetBufferLength ( ) const` `[inline]`

27.211.3.5 `bool gdcm::PNMCodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` `[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.211.3.6 `bool gdcm::PNMCodec::Read ( const char * filename, DataElement & out ) const`

27.211.3.7 `void gdcm::PNMCodec::SetBufferLength ( unsigned long l )` `[inline]`

27.211.3.8 `bool gdcm::PNMCodec::Write ( const char * filename, const DataElement & out ) const`

Examples:

[ExtractIconFromFile.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmPNMCodec.h](#)

## 27.212 gdcm::Preamble Class Reference

DICOM [Preamble](#) (Part 10)

```
#include <gdcmPreamble.h>
```

### Public Member Functions

- [Preamble](#) ()
- [Preamble](#) ([Preamble](#) const &)
- [~Preamble](#) ()
- void [Clear](#) ()
- void [Create](#) ()
- const char \* [GetInternal](#) () const
- [VL GetLength](#) () const
- bool [IsEmpty](#) () const
- [Preamble](#) & [operator=](#) ([Preamble](#) const &)
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [Remove](#) ()
- void [Valid](#) ()
- std::ostream const & [Write](#) (std::ostream &os) const

### Protected Member Functions

- bool [IsValid](#) () const

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Preamble](#) &\_val)

### 27.212.1 Detailed Description

DICOM [Preamble](#) (Part 10)

### 27.212.2 Constructor & Destructor Documentation

27.212.2.1 `gdcM::Preamble::Preamble ( )`

27.212.2.2 `gdcM::Preamble::~~Preamble ( )`

27.212.2.3 `gdcM::Preamble::Preamble ( Preamble const & )` `[inline]`

### 27.212.3 Member Function Documentation

27.212.3.1 `void gdcM::Preamble::Clear ( )`

27.212.3.2 `void gdcM::Preamble::Create ( )`

27.212.3.3 `const char* gdcM::Preamble::GetInternal ( ) const` `[inline]`

27.212.3.4 `VL gdcM::Preamble::GetLength ( ) const` `[inline]`

27.212.3.5 `bool gdcM::Preamble::IsEmpty ( ) const` `[inline]`

27.212.3.6 `bool gdcM::Preamble::IsValid ( ) const` `[inline]`, `[protected]`

27.212.3.7 `Preamble& gdcM::Preamble::operator= ( Preamble const & )` `[inline]`

27.212.3.8 `void gdcM::Preamble::Print ( std::ostream & os ) const`

27.212.3.9 `std::istream& gdcM::Preamble::Read ( std::istream & is )`

27.212.3.10 `void gdcM::Preamble::Remove ( )`

27.212.3.11 `void gdcM::Preamble::Valid ( )`

27.212.3.12 `std::ostream const& gdcM::Preamble::Write ( std::ostream & os ) const`

### 27.212.4 Friends And Related Function Documentation

27.212.4.1 `std::ostream& operator<< ( std::ostream & _os, const Preamble & _val )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcM\\_Preamble.h](#)

## 27.213 gdcM::PresentationContext Class Reference

[PresentationContext](#).

```
#include <gdcM_PresentationContext.h>
```

### Public Types

- `typedef TransferSyntaxArrayType::size_type` [SizeType](#)

- typedef std::vector< std::string > [TransferSyntaxArrayType](#)

## Public Member Functions

- [PresentationContext](#) ()
- [PresentationContext](#) (UIDs::TSName asname, UIDs::TSName tsname=UIDs::ImplicitVRLittleEndianDefault←  
[TransferSyntaxforDICOM](#))
- void [AddTransferSyntax](#) (const char \*tsstr)
- const char \* [GetAbstractSyntax](#) () const
- [SizeType](#) [GetNumberOfTransferSyntaxes](#) () const
- uint8\_t [GetPresentationContextID](#) () const
- const char \* [GetTransferSyntax](#) ([SizeType](#) i) const
- bool [operator==](#) (const [PresentationContext](#) &pc) const
- void [Print](#) (std::ostream &os) const
- void [SetAbstractSyntax](#) (const char \*as)
- void [SetPresentationContextID](#) (uint8\_t id)

## 27.213.1 Detailed Description

[PresentationContext](#).

See also

[PresentationContextAC](#) [PresentationContextRQ](#)

## 27.213.2 Member Typedef Documentation

27.213.2.1 typedef [TransferSyntaxArrayType::size\\_type](#) [gdcm::PresentationContext::SizeType](#)

27.213.2.2 typedef std::vector< std::string> [gdcm::PresentationContext::TransferSyntaxArrayType](#)

## 27.213.3 Constructor & Destructor Documentation

27.213.3.1 [gdcm::PresentationContext::PresentationContext](#) ( )

27.213.3.2 [gdcm::PresentationContext::PresentationContext](#) ( UIDs::TSName asname, UIDs::TSName tsname =  
UIDs::ImplicitVRLittleEndianDefaultTransferSyntaxforDICOM )

Initialize Presentation Context with AbstractSyntax set to asname and with a single [TransferSyntax](#) set to tsname (default to Implicit [VR](#) LittleEndian when not specified ).

## 27.213.4 Member Function Documentation

27.213.4.1 void [gdcm::PresentationContext::AddTransferSyntax](#) ( const char \* tsstr )

27.213.4.2 const char\* [gdcm::PresentationContext::GetAbstractSyntax](#) ( ) const [\[inline\]](#)

27.213.4.3 [SizeType](#) [gdcm::PresentationContext::GetNumberOfTransferSyntaxes](#) ( ) const [\[inline\]](#)

- 27.213.4.4 `uint8_t gdcm::PresentationContext::GetPresentationContextID ( ) const`
- 27.213.4.5 `const char* gdcm::PresentationContext::GetTransferSyntax ( SizeType i ) const` `[inline]`
- 27.213.4.6 `bool gdcm::PresentationContext::operator== ( const PresentationContext & pc ) const` `[inline]`
- 27.213.4.7 `void gdcm::PresentationContext::Print ( std::ostream & os ) const`
- 27.213.4.8 `void gdcm::PresentationContext::SetAbstractSyntax ( const char * as )` `[inline]`
- 27.213.4.9 `void gdcm::PresentationContext::SetPresentationContextID ( uint8_t id )`

The documentation for this class was generated from the following file:

- [gdcmPresentationContext.h](#)

## 27.214 gdcm::network::PresentationContextAC Class Reference

[PresentationContextAC Table](#) 9-18 PRESENTATION CONTEXT ITEM FIELDS.

```
#include <gdcmPresentationContextAC.h>
```

### Public Member Functions

- [PresentationContextAC](#) ()
- `uint8_t GetPresentationContextID () const`
- `uint8_t GetReason () const`
- `TransferSyntaxSub const & GetTransferSyntax () const`
- `void Print (std::ostream &os) const`
- `std::istream & Read (std::istream &is)`
- `void SetPresentationContextID (uint8_t id)`
- `void SetReason (uint8_t r)`
- `void SetTransferSyntax (TransferSyntaxSub const &ts)`
- `size_t Size () const`
- `const std::ostream & Write (std::ostream &os) const`

### 27.214.1 Detailed Description

[PresentationContextAC Table](#) 9-18 PRESENTATION CONTEXT ITEM FIELDS.

See also

[PresentationContext](#)

### 27.214.2 Constructor & Destructor Documentation

- 27.214.2.1 `gdcm::network::PresentationContextAC::PresentationContextAC ( )`

### 27.214.3 Member Function Documentation

- 27.214.3.1 `uint8_t gdcm::network::PresentationContextAC::GetPresentationContextID ( ) const` `[inline]`
- 27.214.3.2 `uint8_t gdcm::network::PresentationContextAC::GetReason ( ) const` `[inline]`
- 27.214.3.3 `TransferSyntaxSub const& gdcm::network::PresentationContextAC::GetTransferSyntax ( ) const` `[inline]`
- 27.214.3.4 `void gdcm::network::PresentationContextAC::Print ( std::ostream & os ) const`
- 27.214.3.5 `std::istream& gdcm::network::PresentationContextAC::Read ( std::istream & is )`
- 27.214.3.6 `void gdcm::network::PresentationContextAC::SetPresentationContextID ( uint8_t id )`
- 27.214.3.7 `void gdcm::network::PresentationContextAC::SetReason ( uint8_t r )` `[inline]`
- 27.214.3.8 `void gdcm::network::PresentationContextAC::SetTransferSyntax ( TransferSyntaxSub const & ts )`
- 27.214.3.9 `size_t gdcm::network::PresentationContextAC::Size ( ) const`
- 27.214.3.10 `const std::ostream& gdcm::network::PresentationContextAC::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:

- [gdcmPresentationContextAC.h](#)

## 27.215 gdcm::PresentationContextGenerator Class Reference

[PresentationContextGenerator](#) This class is responsible for generating the proper [PresentationContext](#) that will be used in subsequent operation during a DICOM Query/Retrieve association. The step of the association is very sensible as special care need to be taken to explicitly define what instance are going to be send and how they are encoded.

```
#include <gdcmPresentationContextGenerator.h>
```

### Public Types

- `typedef std::vector< PresentationContext > PresentationContextArrayType`
- `typedef PresentationContextArrayType::size_type SizeType`

### Public Member Functions

- [PresentationContextGenerator](#) ()
- `bool GenerateFromFileNames (const Directory::FileNamesType &files)`
- `bool GenerateFromUID (UIDs::TSName asname)`  
*Generate the [PresentationContext](#) array from a UID (eg. VerificationSOPClass)*
- `PresentationContextArrayType const & GetPresentationContexts ()`
- `void SetDefaultTransferSyntax (const TransferSyntax &ts)`  
*Not implemented for now. GDCM internally uses Implicit Little Endian.*
- `void SetMergeModeToAbstractSyntax ()`
- `void SetMergeModeToTransferSyntax ()`

## Protected Member Functions

- bool [AddPresentationContext](#) (const char \*as, const char \*ts)
- const char \* [GetDefaultTransferSyntax](#) () const

### 27.215.1 Detailed Description

[PresentationContextGenerator](#) This class is responsible for generating the proper [PresentationContext](#) that will be used in subsequent operation during a DICOM Query/Retrieve association. The step of the association is very sensible as special care need to be taken to explicitly define what instance are going to be send and how they are encoded.

For example a [PresentationContext](#) will express that negotiation requires that CT [Image](#) Storage are send using JPEG Lossless, while US [Image](#) Storage are sent using RLE Transfer Syntax.

Two very different API are exposed one which will always default to little endian transfer syntax see [GenerateFromUID\(\)](#) This API is used for C-ECHO, C-FIND and C-MOVE (SCU). Another API: [GenerateFromFileNames\(\)](#) is used for C-STORE (SCU) as it will loop over all filenames argument to detect the actual encoding. and therefore find the proper encoding to be used.

Two modes are available. The default mode (SetMergeModeToAbstractSyntax) append [PresentationContext](#) (one AbstractSyntax and one [TransferSyntax](#)), as long a they are different. Eg MR [Image](#) Storage/JPEG2000 and MR [Image](#) Storage/JPEGLossless would be considered different. the other mode SetMergeModeToTransferSyntax merge any new [TransferSyntax](#) to the already existing [PresentationContext](#) in order to re-use the same AbstractSyntax.

See also

[PresentationContext](#)

Examples:

[CStoreQtProgress.cxx](#).

### 27.215.2 Member Typedef Documentation

27.215.2.1 `typedef std::vector<PresentationContext> gdcm::PresentationContextGenerator::PresentationContextArrayType`

27.215.2.2 `typedef PresentationContextArrayType::size_type gdcm::PresentationContextGenerator::SizeType`

### 27.215.3 Constructor & Destructor Documentation

27.215.3.1 `gdcm::PresentationContextGenerator::PresentationContextGenerator ( )`

### 27.215.4 Member Function Documentation

27.215.4.1 `bool gdcm::PresentationContextGenerator::AddPresentationContext ( const char * as, const char * ts )`  
[protected]

27.215.4.2 `bool gdcm::PresentationContextGenerator::GenerateFromFileNames ( const Directory::FileNamesType & files )`

Generate the [PresentationContext](#) array from a File-Set. [File](#) specified needs to be valid DICOM files. Used for C-STORE operations



Examples:

[CStoreQtProgress.cxx](#).

27.215.4.3 `bool gdcm::PresentationContextGenerator::GenerateFromUID ( UIDs::TSName asname )`

Generate the [PresentationContext](#) array from a UID (eg. [VerificationSOPClass](#))

27.215.4.4 `const char* gdcm::PresentationContextGenerator::GetDefaultTransferSyntax ( ) const` `[protected]`

27.215.4.5 `PresentationContextArrayType const& gdcm::PresentationContextGenerator::GetPresentationContexts ( )`  
`[inline]`

Examples:

[CStoreQtProgress.cxx](#).

27.215.4.6 `void gdcm::PresentationContextGenerator::SetDefaultTransferSyntax ( const TransferSyntax & ts )`

Not implemented for now. GDCM internally uses Implicit Little Endian.

27.215.4.7 `void gdcm::PresentationContextGenerator::SetMergeModeToAbstractSyntax ( )`

27.215.4.8 `void gdcm::PresentationContextGenerator::SetMergeModeToTransferSyntax ( )`

The documentation for this class was generated from the following file:

- [gdcmPresentationContextGenerator.h](#)

## 27.216 gdcm::network::PresentationContextRQ Class Reference

[PresentationContextRQ](#) Table 9-13 PRESENTATION CONTEXT ITEM FIELDS.

```
#include <gdcmPresentationContextRQ.h>
```

### Public Types

- typedef std::vector< [TransferSyntaxSub](#) >::size\_type [SizeType](#)

### Public Member Functions

- [PresentationContextRQ](#) ( )
- [PresentationContextRQ](#) ([UIDs::TSName](#) *asname*, [UIDs::TSName](#) *tsname*=[UIDs::ImplicitVRLittleEndianDefault](#)↔[TransferSyntaxforDICOM](#))
- [PresentationContextRQ](#) (const [PresentationContext](#) &*pc*)
- void [AddTransferSyntax](#) ([TransferSyntaxSub](#) const &*ts*)
- [AbstractSyntax](#) const & [GetAbstractSyntax](#) ( ) const

- [AbstractSyntax](#) & [GetAbstractSyntax](#) ( )
- [SizeType](#) [GetNumberOfTransferSyntaxes](#) ( ) const
- [uint8\\_t](#) [GetPresentationContextID](#) ( ) const
- [TransferSyntaxSub](#) const & [GetTransferSyntax](#) ( [SizeType](#) i ) const
- [TransferSyntaxSub](#) & [GetTransferSyntax](#) ( [SizeType](#) i )
- [std::vector](#)< [TransferSyntaxSub](#) > const & [GetTransferSyntaxes](#) ( ) const
- [bool](#) [operator==](#) ( const [PresentationContextRQ](#) &pc ) const
- [void](#) [Print](#) ( [std::ostream](#) &os ) const
- [std::istream](#) & [Read](#) ( [std::istream](#) &is )
- [void](#) [SetAbstractSyntax](#) ( [AbstractSyntax](#) const &as )
- [void](#) [SetPresentationContextID](#) ( [uint8\\_t](#) id )
- [size\\_t](#) [Size](#) ( ) const
- [const std::ostream](#) & [Write](#) ( [std::ostream](#) &os ) const

### 27.216.1 Detailed Description

[PresentationContextRQ](#) Table 9-13 PRESENTATION CONTEXT ITEM FIELDS.

See also

[PresentationContextAC](#)

### 27.216.2 Member Typedef Documentation

27.216.2.1 `typedef std::vector<TransferSyntaxSub>::size_type gdcm::network::PresentationContextRQ::SizeType`

### 27.216.3 Constructor & Destructor Documentation

27.216.3.1 `gdcm::network::PresentationContextRQ::PresentationContextRQ ( )`

27.216.3.2 `gdcm::network::PresentationContextRQ::PresentationContextRQ ( UIDs::TSName asname, UIDs::TSName tsname = UIDs::ImplicitVRLittleEndianDefaultTransferSyntaxforDICOM )`

Initialize Presentation Context with [AbstractSyntax](#) set to asname and with a single [TransferSyntax](#) set to tsname (default to Implicit [VR](#) LittleEndian when not specified ).

27.216.3.3 `gdcm::network::PresentationContextRQ::PresentationContextRQ ( const PresentationContext & pc )`

### 27.216.4 Member Function Documentation

27.216.4.1 `void gdcm::network::PresentationContextRQ::AddTransferSyntax ( TransferSyntaxSub const & ts )`

27.216.4.2 `AbstractSyntax const& gdcm::network::PresentationContextRQ::GetAbstractSyntax ( ) const` [[inline](#)]

27.216.4.3 `AbstractSyntax& gdcm::network::PresentationContextRQ::GetAbstractSyntax ( )` [[inline](#)]

27.216.4.4 `SizeType gdcm::network::PresentationContextRQ::GetNumberOfTransferSyntaxes ( ) const` [[inline](#)]

27.216.4.5 `uint8\_t gdcm::network::PresentationContextRQ::GetPresentationContextID ( ) const`

- 27.216.4.6 **TransferSyntaxSub** const& gdcm::network::PresentationContextRQ::GetTransferSyntax ( **SizeType** *i* ) const  
[inline]
- 27.216.4.7 **TransferSyntaxSub**& gdcm::network::PresentationContextRQ::GetTransferSyntax ( **SizeType** *i* ) [inline]
- 27.216.4.8 **std::vector<TransferSyntaxSub>** const& gdcm::network::PresentationContextRQ::GetTransferSyntaxes ( ) const  
[inline]
- 27.216.4.9 **bool** gdcm::network::PresentationContextRQ::operator== ( const **PresentationContextRQ** & *pc* ) const  
[inline]
- 27.216.4.10 **void** gdcm::network::PresentationContextRQ::Print ( **std::ostream** & *os* ) const
- 27.216.4.11 **std::istream**& gdcm::network::PresentationContextRQ::Read ( **std::istream** & *is* )
- 27.216.4.12 **void** gdcm::network::PresentationContextRQ::SetAbstractSyntax ( **AbstractSyntax** const & *as* )
- 27.216.4.13 **void** gdcm::network::PresentationContextRQ::SetPresentationContextID ( **uint8\_t** *id* )
- 27.216.4.14 **size\_t** gdcm::network::PresentationContextRQ::Size ( ) const
- 27.216.4.15 **const std::ostream**& gdcm::network::PresentationContextRQ::Write ( **std::ostream** & *os* ) const

The documentation for this class was generated from the following file:

- [gdcmPresentationContextRQ.h](#)

## 27.217 gdcm::network::PresentationDataValue Class Reference

[PresentationDataValue](#) Table 9-23 PRESENTATION-DATA-VALUE ITEM FIELDS.

```
#include <gdcmPresentationDataValue.h>
```

### Public Member Functions

- [PresentationDataValue](#) ()
- **const std::string** & [GetBlob](#) () const
- **bool** [GetIsCommand](#) () const
- **bool** [GetIsLastFragment](#) () const
- **uint8\_t** [GetMessageHeader](#) () const
- **uint8\_t** [GetPresentationContextID](#) () const
- **void** [Print](#) (**std::ostream** & *os*) const
- **std::istream** & [Read](#) (**std::istream** & *is*)
- **std::istream** & [ReadInto](#) (**std::istream** & *is*, **std::ostream** & *os*)
- **void** [SetBlob](#) (**const std::string** & *partialblob*)
- **void** [SetCommand](#) (**bool** *inCommand*)
- **void** [SetDataSet](#) (**const DataSet** & *ds*)
- **void** [SetLastFragment](#) (**bool** *inLast*)
- **void** [SetMessageHeader](#) (**uint8\_t** *messageheader*)
- **void** [SetPresentationContextID](#) (**uint8\_t** *id*)
- **size\_t** [Size](#) () const
- **const std::ostream** & [Write](#) (**std::ostream** & *os*) const

## Static Public Member Functions

- static [DataSet ConcatenatePDVBlobs](#) (const std::vector< [PresentationDataValue](#) > &inPDVs)
- static [DataSet ConcatenatePDVBlobsAsExplicit](#) (const std::vector< [PresentationDataValue](#) > &inPDVs)

### 27.217.1 Detailed Description

[PresentationDataValue Table](#) 9-23 PRESENTATION-DATA-VALUE ITEM FIELDS.

### 27.217.2 Constructor & Destructor Documentation

27.217.2.1 `gdcm::network::PresentationDataValue::PresentationDataValue ( )`

### 27.217.3 Member Function Documentation

27.217.3.1 `static DataSet gdcm::network::PresentationDataValue::ConcatenatePDVBlobs ( const std::vector< PresentationDataValue > & inPDVs ) [static]`

Warning

[DataSet](#) will be read as Implicit Little Endian TS

27.217.3.2 `static DataSet gdcm::network::PresentationDataValue::ConcatenatePDVBlobsAsExplicit ( const std::vector< PresentationDataValue > & inPDVs ) [static]`

27.217.3.3 `const std::string& gdcm::network::PresentationDataValue::GetBlob ( ) const`

27.217.3.4 `bool gdcm::network::PresentationDataValue::GetIsCommand ( ) const`

27.217.3.5 `bool gdcm::network::PresentationDataValue::GetIsLastFragment ( ) const`

27.217.3.6 `uint8_t gdcm::network::PresentationDataValue::GetMessageHeader ( ) const [inline]`

27.217.3.7 `uint8_t gdcm::network::PresentationDataValue::GetPresentationContextID ( ) const [inline]`

27.217.3.8 `void gdcm::network::PresentationDataValue::Print ( std::ostream & os ) const`

27.217.3.9 `std::istream& gdcm::network::PresentationDataValue::Read ( std::istream & is )`

27.217.3.10 `std::istream& gdcm::network::PresentationDataValue::ReadInto ( std::istream & is, std::ostream & os )`

27.217.3.11 `void gdcm::network::PresentationDataValue::SetBlob ( const std::string & partialblob )`

27.217.3.12 `void gdcm::network::PresentationDataValue::SetCommand ( bool inCommand )`

27.217.3.13 `void gdcm::network::PresentationDataValue::SetDataSet ( const DataSet & ds )`

Set [DataSet](#). Write [DataSet](#) in implicit.

## Warning

size of dataset should be below maxpdusize

27.217.3.14 void gdcm::network::PresentationDataValue::SetLastFragment ( bool *inLast* )

27.217.3.15 void gdcm::network::PresentationDataValue::SetMessageHeader ( uint8\_t *messageheader* ) [inline]

27.217.3.16 void gdcm::network::PresentationDataValue::SetPresentationContextID ( uint8\_t *id* ) [inline]

27.217.3.17 size\_t gdcm::network::PresentationDataValue::Size ( ) const

27.217.3.18 const std::ostream& gdcm::network::PresentationDataValue::Write ( std::ostream & *os* ) const

The documentation for this class was generated from the following file:


- [gdcmPresentationDataValue.h](#)

## 27.218 gdcm::Printer Class Reference


[Printer](#) class.

```
#include <gdcmPrinter.h>
```

Inheritance diagram for gdcm::Printer:

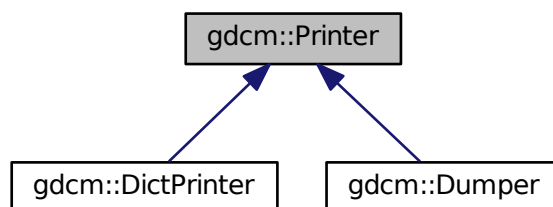

Collaboration diagram for `gdcM::Printer`:

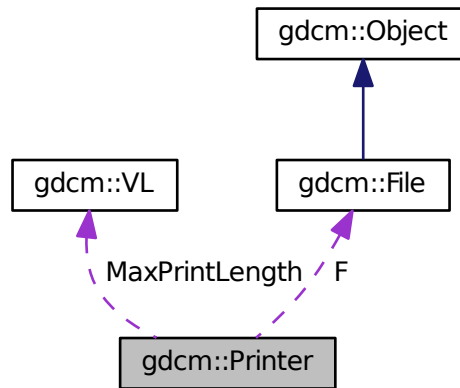

## Public Types

- enum `PrintStyles` {  
`VERBOSE_STYLE` = 0,  
`CONDENSED_STYLE`,  
`XML` }

## Public Member Functions

- `Printer ()`
- `~Printer ()`
- `PrintStyles GetPrintStyle ()` const  
*Get PrintStyle value.*
- `void Print (std::ostream &os)`  
*Print.*
- `void PrintDataSet (const DataSet &ds, std::ostream &os, const std::string &s="")`  
*Print an individual dataset.*
- `void SetColor (bool c)`  
*Set color mode or not.*
- `void SetFile (File const &f)`  
*Set file.*
- `void SetStyle (PrintStyles ps)`  
*Set PrintStyle value.*

## Protected Member Functions

- `VR PrintDataElement (std::ostringstream &os, const Dicts &dicts, const DataSet &ds, const DataElement &de, std::ostream &out, std::string const &indent)`

- void [PrintSQ](#) (const [SequenceOfItems](#) \*sqi, std::ostream &os, std::string const &indent)

## Protected Attributes

- const [File](#) \* F
- [VL MaxPrintLength](#)
- [PrintStyles](#) [PrintStyle](#)

## 27.218.1 Detailed Description

[Printer](#) class.

## 27.218.2 Member Enumeration Documentation

### 27.218.2.1 enum gdcm::Printer::PrintStyles

Enumerator

***VERBOSE\_STYLE***

***CONDENSED\_STYLE***

***XML***

## 27.218.3 Constructor & Destructor Documentation

### 27.218.3.1 gdcm::Printer::Printer ( )

### 27.218.3.2 gdcm::Printer::~~Printer ( )

## 27.218.4 Member Function Documentation

### 27.218.4.1 [PrintStyles](#) gdcm::Printer::GetPrintStyle ( ) const [inline]

Get PrintStyle value.

### 27.218.4.2 void gdcm::Printer::Print ( std::ostream & os )

Print.

### 27.218.4.3 VR gdcm::Printer::PrintDataElement ( std::ostream & os, const [Dicts](#) & *dicts*, const [DataSet](#) & *ds*, const [DataElement](#) & *de*, std::ostream & *out*, std::string const & *indent* ) [protected]

### 27.218.4.4 void gdcm::Printer::PrintDataSet ( const [DataSet](#) & *ds*, std::ostream & *os*, const std::string & *s* = " " )

Print an individual dataset.

27.218.4.5 void `gdcmm::Printer::PrintSQ` ( const `SequenceOfItems` \* *sqi*, `std::ostream` & *os*, `std::string` const & *indent* )  
[protected]

27.218.4.6 void `gdcmm::Printer::SetColor` ( bool *c* )

Set color mode or not.

27.218.4.7 void `gdcmm::Printer::SetFile` ( `File` const & *f* ) [inline]

Set file.

27.218.4.8 void `gdcmm::Printer::SetStyle` ( `PrintStyles` *ps* ) [inline]

Set `PrintStyle` value.

## 27.218.5 Member Data Documentation

27.218.5.1 const `File`\* `gdcmm::Printer::F` [protected]

27.218.5.2 VL `gdcmm::Printer::MaxPrintLength` [protected]

27.218.5.3 `PrintStyles` `gdcmm::Printer::PrintStyle` [protected]

The documentation for this class was generated from the following file:

- [gdcmmPrinter.h](#)

## 27.219 gdcmm::PrivateDict Class Reference

Private [Dict](#).

```
#include <gdcmmDict.h>
```

### Public Member Functions

- [PrivateDict](#) ()
- [~PrivateDict](#) ()
- void [AddDictEntry](#) (const [PrivateTag](#) &tag, const [DictEntry](#) &de)
- bool [FindDictEntry](#) (const [PrivateTag](#) &tag) const
- const [DictEntry](#) & [GetDictEntry](#) (const [PrivateTag](#) &tag) const
- bool [IsEmpty](#) () const
- void [PrintXML](#) () const
- bool [RemoveDictEntry](#) (const [PrivateTag](#) &tag)

### Protected Member Functions

- void [LoadDefault](#) ()



## Friends

- class [Dicts](#)
- `std::ostream & operator<< (std::ostream &os, const PrivateDict &val)`

## 27.219.1 Detailed Description

Private [Dict](#).

## 27.219.2 Constructor & Destructor Documentation

27.219.2.1 `gdcm::PrivateDict::PrivateDict ( )` `[inline]`

27.219.2.2 `gdcm::PrivateDict::~~PrivateDict ( )` `[inline]`

## 27.219.3 Member Function Documentation

27.219.3.1 `void gdcm::PrivateDict::AddDictEntry ( const PrivateTag &tag, const DictEntry &de )` `[inline]`

References `gdcm::DictEntry::GetVM()`, `gdcm::DictEntry::GetVR()`, `gdcm::DictEntry::SetVR()`, and `gdcm::VR::UN`.

27.219.3.2 `bool gdcm::PrivateDict::FindDictEntry ( const PrivateTag &tag ) const` `[inline]`

27.219.3.3 `const DictEntry& gdcm::PrivateDict::GetDictEntry ( const PrivateTag &tag ) const` `[inline]`

27.219.3.4 `bool gdcm::PrivateDict::IsEmpty ( ) const` `[inline]`

27.219.3.5 `void gdcm::PrivateDict::LoadDefault ( )` `[protected]`

27.219.3.6 `void gdcm::PrivateDict::PrintXML ( ) const` `[inline]`

References `gdcm::Tag::GetElement()`, `gdcm::Tag::GetGroup()`, `gdcm::DictEntry::GetName()`, `gdcm::PrivateTag::GetOwner()`, `gdcm::DictEntry::GetVM()`, and `gdcm::DictEntry::GetVR()`.

27.219.3.7 `bool gdcm::PrivateDict::RemoveDictEntry ( const PrivateTag &tag )` `[inline]`

Remove entry 'tag'. Return true on success (element was found and remove). return false if element was not found.

## 27.219.4 Friends And Related Function Documentation

27.219.4.1 `friend class Dicts` `[friend]`

27.219.4.2 `std::ostream& operator<< ( std::ostream &os, const PrivateDict &val )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmDict.h](#)

## 27.220 gdcmm::PrivateTag Class Reference


Class to represent a Private DICOM Data [Element](#) ([Attribute](#)) [Tag](#) (Group, [Element](#), Owner)

```
#include <gdcmmPrivateTag.h>
```

Inheritance diagram for gdcmm::PrivateTag:

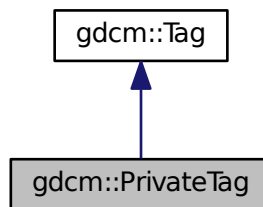

Collaboration diagram for gdcmm::PrivateTag:

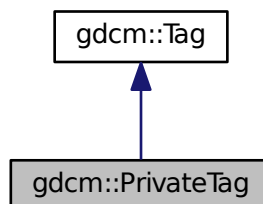

### Public Member Functions

- [PrivateTag](#) (uint16\_t group=0, uint16\_t element=0, const char \*owner="")
- [PrivateTag](#) ([Tag](#) const &t, const char \*owner="")
- [DataElement](#) [GetAsDataElement](#) () const
- const char \* [GetOwner](#) () const
- bool [operator<](#) (const [PrivateTag](#) &\_val) const
- bool [ReadFromCommaSeparatedString](#) (const char \*str)
- void [SetOwner](#) (const char \*owner)

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [PrivateTag](#) &\_val)

### 27.220.1 Detailed Description

Class to represent a Private DICOM Data [Element](#) ([Attribute](#)) [Tag](#) (Group, [Element](#), Owner)

#### Note

private tag have element value in: [0x10,0xff], for instance 0x0009,0x0000 is NOT a private tag

#### Examples:

[ChangePrivateTags.cxx](#), [csa2img.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [ELSCINT1WaveToText.cxx](#), [GetSubSequenceData.cxx](#), [i↔U22tomultisc.cxx](#), [MrProtocol.cxx](#), [pmsct\\_rgb1.cxx](#), [PublicDict.cxx](#), [ReadGEMSSDO.cxx](#), and [rle2img.cxx](#).

### 27.220.2 Constructor & Destructor Documentation

27.220.2.1 `gdcm::PrivateTag::PrivateTag ( uint16_t group = 0, uint16_t element = 0, const char * owner = " " )` `[inline]`

27.220.2.2 `gdcm::PrivateTag::PrivateTag ( Tag const & t, const char * owner = " " )` `[inline]`

References `gdcm::Tag::GetElement()`.

### 27.220.3 Member Function Documentation

27.220.3.1 `DataElement gdcm::PrivateTag::GetAsDataElement ( )` `const`

27.220.3.2 `const char* gdcm::PrivateTag::GetOwner ( )` `const` `[inline]`

#### Examples:

[PublicDict.cxx](#).

Referenced by `gdcm::PrivateDict::PrintXML()`.

27.220.3.3 `bool gdcm::PrivateTag::operator< ( const PrivateTag & _val )` `const`

27.220.3.4 `bool gdcm::PrivateTag::ReadFromCommaSeparatedString ( const char * str )`

Read [PrivateTag](#) from a string. [Element](#) number will be truncated to 8bits. Eg: "1234,5678,GDCM" is private tag: (1234,78,"GDCM")

27.220.3.5 `void gdcm::PrivateTag::SetOwner ( const char * owner )` `[inline]`

### 27.220.4 Friends And Related Function Documentation

27.220.4.1 `std::ostream& operator<< ( std::ostream & _os, const PrivateTag & _val )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmPrivateTag.h](#)

## 27.221 gdcM::ProgressEvent Class Reference


[ProgressEvent](#) Special type of event triggered during.

```
#include <gdcMProgressEvent.h>
```

Inheritance diagram for gdcM::ProgressEvent:

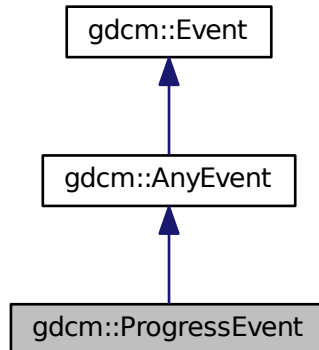

Collaboration diagram for gdcM::ProgressEvent:

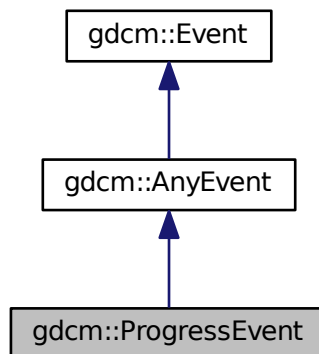

### Public Types

- typedef [ProgressEvent](#) Self
- typedef [AnyEvent](#) Superclass

## Public Member Functions

- [ProgressEvent](#) (double p=0)
- [ProgressEvent](#) (const [Self](#) &s)
- virtual [~ProgressEvent](#) ()
- virtual bool [CheckEvent](#) (const [::gdcm::Event](#) \*e) const
- virtual const char \* [GetEventName](#) () const
- double [GetProgress](#) () const
- virtual [::gdcm::Event](#) \* [MakeObject](#) () const
- void [SetProgress](#) (double p)

### 27.221.1 Detailed Description

[ProgressEvent](#) Special type of event triggered during.

See also

[AnyEvent](#)

### 27.221.2 Member Typedef Documentation

27.221.2.1 typedef [ProgressEvent](#) [gdcm::ProgressEvent::Self](#)

27.221.2.2 typedef [AnyEvent](#) [gdcm::ProgressEvent::Superclass](#)

### 27.221.3 Constructor & Destructor Documentation

27.221.3.1 [gdcm::ProgressEvent::ProgressEvent](#) ( double *p* = 0 ) [\[inline\]](#)

27.221.3.2 virtual [gdcm::ProgressEvent::~~ProgressEvent](#) ( ) [\[inline\]](#),[\[virtual\]](#)

27.221.3.3 [gdcm::ProgressEvent::ProgressEvent](#) ( const [Self](#) & *s* ) [\[inline\]](#)

### 27.221.4 Member Function Documentation

27.221.4.1 virtual bool [gdcm::ProgressEvent::CheckEvent](#) ( const [::gdcm::Event](#) \* *e* ) const [\[inline\]](#),[\[virtual\]](#)

27.221.4.2 virtual const char\* [gdcm::ProgressEvent::GetEventName](#) ( ) const [\[inline\]](#),[\[virtual\]](#)

Return the StringName associated with the event.

Implements [gdcm::Event](#).

27.221.4.3 double [gdcm::ProgressEvent::GetProgress](#) ( ) const [\[inline\]](#)

27.221.4.4 virtual [::gdcm::Event](#)\* [gdcm::ProgressEvent::MakeObject](#) ( ) const [\[inline\]](#),[\[virtual\]](#)

Create an [Event](#) of this type This method work as a Factory for creating events of each particular type.

Implements [gdcm::Event](#).

27.221.4.5 void `gdcm::ProgressEvent::SetProgress ( double p )` [inline]

The documentation for this class was generated from the following file:


- [gdcmProgressEvent.h](#)

## 27.222 `gdcm::PVRGCodec` Class Reference


[PVRGCodec](#).

```
#include <gdcmPVRGCodec.h>
```

Inheritance diagram for `gdcm::PVRGCodec`:

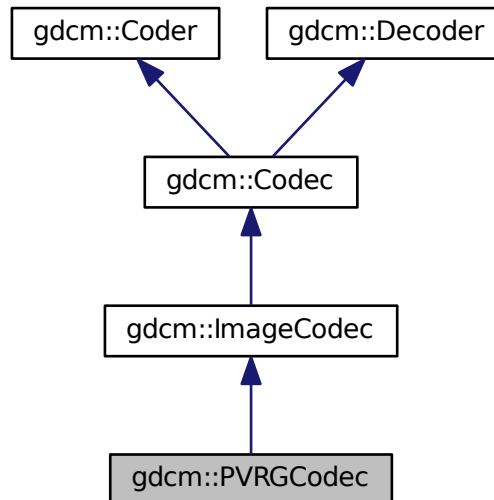

Collaboration diagram for `gdcm::PVRGCodec`:

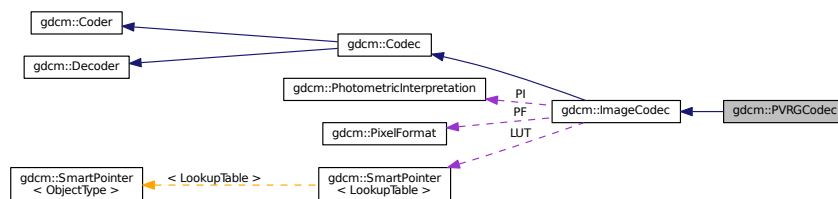

## Public Member Functions

- [PVRGCodec](#) ()
- [~PVRGCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &ts) const  
*Return whether this coder support this transfer syntax (can code it)*
- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const  
*Return whether this decoder support this transfer syntax (can decode it)*
- virtual [ImageCodec](#) \* [Clone](#) () const
- bool [Code](#) ([DataElement](#) const &in, [DataElement](#) &out)  
*Code.*
- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)  
*Decode.*
- void [SetLossyFlag](#) (bool l)

## Additional Inherited Members

### 27.222.1 Detailed Description

[PVRGCodec](#).

#### Note

pvr is a broken implementation of the JPEG standard. It is known to have a bug in the 16bits lossless implementation of the standard.

In an ideal world, you should not need this codec at all. But to support some broken file such as:

PHILIPS\_Gyroscan-12-Jpeg\_Extended\_Process\_2\_4.dcm

we have to...

### 27.222.2 Constructor & Destructor Documentation

27.222.2.1 `gdcm::PVRGCodec::PVRGCodec ( )`

27.222.2.2 `gdcm::PVRGCodec::~~PVRGCodec ( )`

### 27.222.3 Member Function Documentation

27.222.3.1 `bool gdcm::PVRGCodec::CanCode ( TransferSyntax const & ) const` [virtual]

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.222.3.2 `bool gdcm::PVRGCodec::CanDecode ( TransferSyntax const & ) const` [virtual]

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.222.3.3 `virtual ImageCodec* gdcm::PVRGCodec::Clone ( ) const` [virtual]

Implements [gdcm::ImageCodec](#).

27.222.3.4 `bool gdcm::PVRGCodec::Code ( DataElement const & in_, DataElement & out_ )` [virtual]

Code.

Reimplemented from [gdcm::Coder](#).

27.222.3.5 `bool gdcm::PVRGCodec::Decode ( DataElement const &, DataElement & )` [virtual]

Decode.

Reimplemented from [gdcm::ImageCodec](#).

27.222.3.6 `void gdcm::PVRGCodec::SetLossyFlag ( bool l )`

The documentation for this class was generated from the following file:

- [gdcmPVRGCodec.h](#)

## 27.223 gdcm::PythonFilter Class Reference

[PythonFilter](#) [PythonFilter](#) is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.

```
#include <gdcmPythonFilter.h>
```

### Public Member Functions

- [PythonFilter](#) ()
- [~PythonFilter](#) ()
- [File](#) & [GetFile](#) ()
- const [File](#) & [GetFile](#) () const
- void [SetDicts](#) (const [Dicts](#) &dicts)
- void [SetFile](#) (const [File](#) &f)
- PyObject \* [ToPyObject](#) (const [Tag](#) &t) const
- void [UseDictAlways](#) (bool use)

### 27.223.1 Detailed Description

[PythonFilter](#) [PythonFilter](#) is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.



## 27.223.2 Constructor & Destructor Documentation

27.223.2.1 `gdcm::PythonFilter::PythonFilter ( )`

27.223.2.2 `gdcm::PythonFilter::~~PythonFilter ( )`

## 27.223.3 Member Function Documentation

27.223.3.1 `File& gdcm::PythonFilter::GetFile ( )` `[inline]`

27.223.3.2 `const File& gdcm::PythonFilter::GetFile ( ) const` `[inline]`

27.223.3.3 `void gdcm::PythonFilter::SetDicts ( const Dicts & dicts )`

27.223.3.4 `void gdcm::PythonFilter::SetFile ( const File & f )` `[inline]`

27.223.3.5 `PyObject* gdcm::PythonFilter::ToPyObject ( const Tag & t ) const`

27.223.3.6 `void gdcm::PythonFilter::UseDictAlways ( bool use )` `[inline]`

The documentation for this class was generated from the following file:


- [gdcmPythonFilter.h](#)

## 27.224 gdcm::QueryBase Class Reference

[QueryBase](#) contains: the base class for constructing a query dataset for a C-FIND and a C-MOVE.

```
#include <gdcmQueryBase.h>
```

Inheritance diagram for `gdcm::QueryBase`:

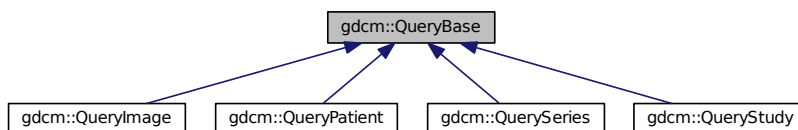

## Public Member Functions

- virtual `~QueryBase ( )`
- `std::vector< Tag > GetAllRequiredTags (const ERootType &inRootType) const`
- `std::vector< Tag > GetAllTags (const ERootType &inRootType) const`
- virtual `std::vector< Tag > GetHierarchicalSearchTags (const ERootType &inRootType) const =0`  
*Return all Unique Key for a particular Query Root type (from the same level and above).*
- virtual `const char * GetName ( ) const =0`
- virtual `std::vector< Tag > GetOptionalTags (const ERootType &inRootType) const =0`

- virtual [DataElement GetQueryLevel](#) () const =0
- virtual std::vector< [Tag](#) > [GetRequiredTags](#) (const [ERootType](#) &inRootType) const =0
- virtual std::vector< [Tag](#) > [GetUniqueTags](#) (const [ERootType](#) &inRootType) const =0

### 27.224.1 Detailed Description

[QueryBase](#) contains: the base class for constructing a query dataset for a C-FIND and a C-MOVE.

There are four levels of C-FIND and C-MOVE query:

- [Patient](#)
- [Study](#)
- [Series](#)
- [Image](#)

Each one has its own required and optional tags. This class provides an interface for getting those tags. This is an interface class.

See 3.4 C 6.1 and 3.4 C 6.2 for the patient and study root query types. These sections define the tags allowed by a particular query. The caller must pass in which root type they want, patient or study. A third root type, Modality Worklist Query, isn't yet supported.

This class (or rather it's derived classes) will be held in the RootQuery types. These query types actually make the dataset, and will use this dataset to list the required, unique, and optional tags for each type of query. This design is somewhat overly complicated, but is kept so that if we ever wanted to try to guess the query type from the given tags, we could do so.

### 27.224.2 Constructor & Destructor Documentation

27.224.2.1 virtual [gdcmm::QueryBase::~~QueryBase](#) ( ) [inline],[virtual]

### 27.224.3 Member Function Documentation

27.224.3.1 std::vector<[Tag](#)> [gdcmm::QueryBase::GetAllRequiredTags](#) ( const [ERootType](#) & *inRootType* ) const

In order to validate a query dataset we need to check that there exists at least one required (or unique) key

27.224.3.2 std::vector<[Tag](#)> [gdcmm::QueryBase::GetAllTags](#) ( const [ERootType](#) & *inRootType* ) const

In order to validate a query dataset, just check for the presence of a tag, not it's requirement level in the spec

27.224.3.3 virtual std::vector<[Tag](#)> [gdcmm::QueryBase::GetHierarchicalSearchTags](#) ( const [ERootType](#) & *inRootType* ) const  
[pure virtual]

Return all Unique Key for a particular Query Root type (from the same level and above).

Implemented in [gdcmm::QueryImage](#), [gdcmm::QueryPatient](#), [gdcmm::QuerySeries](#), and [gdcmm::QueryStudy](#).

27.224.3.4 `virtual const char* gdcm::QueryBase::GetName ( ) const` [pure virtual]

Implemented in [gdcm::QueryImage](#), [gdcm::QueryPatient](#), [gdcm::QuerySeries](#), and [gdcm::QueryStudy](#).

27.224.3.5 `virtual std::vector<Tag> gdcm::QueryBase::GetOptionalTags ( const ERootType & inRootType ) const` [pure virtual]

Implemented in [gdcm::QueryImage](#), [gdcm::QueryPatient](#), [gdcm::QuerySeries](#), and [gdcm::QueryStudy](#).

27.224.3.6 `virtual DataElement gdcm::QueryBase::GetQueryLevel ( ) const` [pure virtual]

Implemented in [gdcm::QueryImage](#), [gdcm::QueryPatient](#), [gdcm::QuerySeries](#), and [gdcm::QueryStudy](#).

27.224.3.7 `virtual std::vector<Tag> gdcm::QueryBase::GetRequiredTags ( const ERootType & inRootType ) const` [pure virtual]

Implemented in [gdcm::QueryImage](#), [gdcm::QueryPatient](#), [gdcm::QuerySeries](#), and [gdcm::QueryStudy](#).

27.224.3.8 `virtual std::vector<Tag> gdcm::QueryBase::GetUniqueTags ( const ERootType & inRootType ) const` [pure virtual]

Implemented in [gdcm::QueryImage](#), [gdcm::QueryPatient](#), [gdcm::QuerySeries](#), and [gdcm::QueryStudy](#).

The documentation for this class was generated from the following file:

- [gdcmQueryBase.h](#)

## 27.225 gdcm::QueryFactory Class Reference

QueryFactory.h.

```
#include <gdcmQueryFactory.h>
```

### Static Public Member Functions

- static [ECharSet](#) [GetCharacterFromCurrentLocale](#) ()
- static void [ListCharSets](#) (std::ostream &os)  
*List all possible CharSet.*
- static [DataElement](#) [ProduceCharacterSetDataElement](#) (const std::vector< [ECharSet](#) > &inCharSetType)
- static [BaseRootQuery](#) \* [ProduceQuery](#) ([ERootType](#) inRootType, [EQueryType](#) inQueryType, [EQueryLevel](#) inQueryLevel)

### 27.225.1 Detailed Description

QueryFactory.h.

**Note**

contains: a class to produce a query based off of user-entered information

Essentially, this class is used to construct a query based off of user input (typically from the command line; if in code directly, the query itself could just be instantiated)

In theory, could also be used as the interface to validate incoming datasets as belonging to a particular query style

**27.225.2 Member Function Documentation****27.225.2.1 static ECharSet gdcm::QueryFactory::GetCharacterFromCurrentLocale ( ) [static]**

This function will return the corresponding ECharSet associated with the current locale of the running system (based on the value of locale() ).

**27.225.2.2 static void gdcm::QueryFactory::ListCharSets ( std::ostream & os ) [static]**

List all possible CharSet.

**27.225.2.3 static DataElement gdcm::QueryFactory::ProduceCharacterSetDataElement ( const std::vector< ECharSet > & inCharSetType ) [static]**

This function will produce the appropriate dataelement given a list of charsets. The first charset will be used directly, while the second and subsequent will be prepended with "ISO2022 ". Redundant character sets are not permitted, so if they are encountered, they will just be skipped. if UTF8 or GB18030 is used, no subsequent character sets will be used if the vector passed in is empty, then the dataelement that's passed out will be empty and Latin1 is the presumed encoding

**27.225.2.4 static BaseRootQuery\* gdcm::QueryFactory::ProduceQuery ( ERootType inRootType, EQueryType inQueryType, EQueryLevel inQueryLevel ) [static]**

this function will produce a query (basically, a wrapper to a dataset that can validate whether or not the query is a valid cfind/cmove query) and the level of the query (patient, study, series, image). If the user provides an invalid instantiation (ie, study root type, query level of patient), then the result is NULL.

The documentation for this class was generated from the following file:


- [gdcmQueryFactory.h](#)

**27.226 gdcm::QueryImage Class Reference**


[QueryImage](#) contains: class to construct an image-based query for C-FIND and C-MOVE.

```
#include <gdcmQueryImage.h>
```

Inheritance diagram for gdcm::QueryImage:

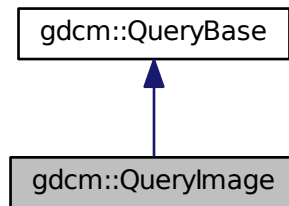

Collaboration diagram for gdcm::QueryImage:

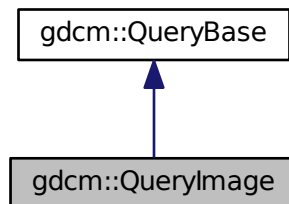

## Public Member Functions

- `std::vector< Tag > GetHierachicalSearchTags (const ERootType &inRootType) const`  
*Return all Unique Key for a particular Query Root type (from the same level and above).*
- `const char * GetName () const`
- `std::vector< Tag > GetOptionalTags (const ERootType &inRootType) const`
- `DataElement GetQueryLevel () const`
- `std::vector< Tag > GetRequiredTags (const ERootType &inRootType) const`
- `std::vector< Tag > GetUniqueTags (const ERootType &inRootType) const`

### 27.226.1 Detailed Description

`QueryImage` contains: class to construct an image-based query for C-FIND and C-MOVE.

### 27.226.2 Member Function Documentation

27.226.2.1 `std::vector<Tag> gdcm::QueryImage::GetHierachicalSearchTags ( const ERootType & inRootType ) const` [virtual]

Return all Unique Key for a particular Query Root type (from the same level and above).

Implements [gdcm::QueryBase](#).

27.226.2.2 `const char* gdcm::QueryImage::GetName ( ) const` [virtual]

Implements [gdcm::QueryBase](#).

27.226.2.3 `std::vector<Tag> gdcm::QueryImage::GetOptionalTags ( const ERootType & inRootType ) const` [virtual]

Implements [gdcm::QueryBase](#).

27.226.2.4 `DataElement gdcm::QueryImage::GetQueryLevel ( ) const` [virtual]

Implements [gdcm::QueryBase](#).

27.226.2.5 `std::vector<Tag> gdcm::QueryImage::GetRequiredTags ( const ERootType & inRootType ) const` [virtual]

Implements [gdcm::QueryBase](#).

27.226.2.6 `std::vector<Tag> gdcm::QueryImage::GetUniqueTags ( const ERootType & inRootType ) const` [virtual]

Implements [gdcm::QueryBase](#).

The documentation for this class was generated from the following file:


- [gdcmQueryImage.h](#)

## 27.227 gdcm::QueryPatient Class Reference


[QueryPatient](#) contains: class to construct a patient-based query for c-find and c-move.

```
#include <gdcmQueryPatient.h>
```

Inheritance diagram for gdcm::QueryPatient:

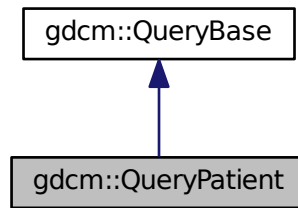

Collaboration diagram for gdcm::QueryPatient:

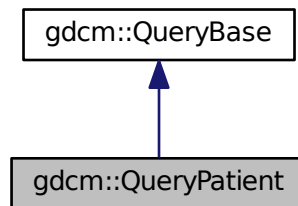

## Public Member Functions

- `std::vector< Tag > GetHierachicalSearchTags (const ERootType &inRootType) const`  
*Return all Unique Key for a particular Query Root type (from the same level and above).*
- `const char * GetName () const`
- `std::vector< Tag > GetOptionalTags (const ERootType &inRootType) const`
- `DataElement GetQueryLevel () const`
- `std::vector< Tag > GetRequiredTags (const ERootType &inRootType) const`
- `std::vector< Tag > GetUniqueTags (const ERootType &inRootType) const`

## 27.227.1 Detailed Description

`QueryPatient` contains: class to construct a patient-based query for c-find and c-move.

## 27.227.2 Member Function Documentation

**27.227.2.1** `std::vector<Tag> gdcm::QueryPatient::GetHierachicalSearchTags ( const ERootType & inRootType ) const` `[virtual]`

Return all Unique Key for a particular Query Root type (from the same level and above).

Implements [gdcm::QueryBase](#).

**27.227.2.2** `const char* gdcm::QueryPatient::GetName ( ) const` `[virtual]`

Implements [gdcm::QueryBase](#).

**27.227.2.3** `std::vector<Tag> gdcm::QueryPatient::GetOptionalTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcm::QueryBase](#).

**27.227.2.4** `DataElement gdcm::QueryPatient::GetQueryLevel ( ) const` `[virtual]`

Implements [gdcm::QueryBase](#).

**27.227.2.5** `std::vector<Tag> gdcm::QueryPatient::GetRequiredTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcm::QueryBase](#).

**27.227.2.6** `std::vector<Tag> gdcm::QueryPatient::GetUniqueTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcm::QueryBase](#).

The documentation for this class was generated from the following file:

- [gdcmQueryPatient.h](#)


## 27.228 gdcm::QuerySeries Class Reference

[QuerySeries](#) contains: class to construct a series-based query for c-find and c-move.


```
#include <gdcmQuerySeries.h>
```



Inheritance diagram for gdcm::QuerySeries:

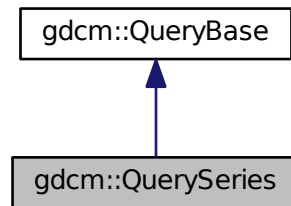

Collaboration diagram for gdcm::QuerySeries:

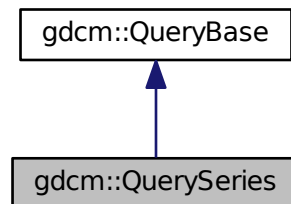

## Public Member Functions

- `std::vector< Tag > GetHierachicalSearchTags (const ERootType &inRootType) const`  
*Return all Unique Key for a particular Query Root type (from the same level and above).*
- `const char * GetName () const`
- `std::vector< Tag > GetOptionalTags (const ERootType &inRootType) const`
- `DataElement GetQueryLevel () const`
- `std::vector< Tag > GetRequiredTags (const ERootType &inRootType) const`
- `std::vector< Tag > GetUniqueTags (const ERootType &inRootType) const`

### 27.228.1 Detailed Description

[QuerySeries](#) contains: class to construct a series-based query for c-find and c-move.

### 27.228.2 Member Function Documentation

27.228.2.1 `std::vector<Tag> gdcM::QuerySeries::GetHierachicalSearchTags ( const ERootType & inRootType ) const` `[virtual]`

Return all Unique Key for a particular Query Root type (from the same level and above).

Implements [gdcM::QueryBase](#).

27.228.2.2 `const char* gdcM::QuerySeries::GetName ( ) const` `[virtual]`

Implements [gdcM::QueryBase](#).

27.228.2.3 `std::vector<Tag> gdcM::QuerySeries::GetOptionalTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcM::QueryBase](#).

27.228.2.4 `DataElement gdcM::QuerySeries::GetQueryLevel ( ) const` `[virtual]`

Implements [gdcM::QueryBase](#).

27.228.2.5 `std::vector<Tag> gdcM::QuerySeries::GetRequiredTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcM::QueryBase](#).

27.228.2.6 `std::vector<Tag> gdcM::QuerySeries::GetUniqueTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcM::QueryBase](#).

The documentation for this class was generated from the following file:


- [gdcMQuerySeries.h](#)

## 27.229 gdcM::QueryStudy Class Reference


QueryStudy.h contains: class to construct a study-based query for C-FIND and C-MOVE.

```
#include <gdcMQueryStudy.h>
```

Inheritance diagram for gdcm::QueryStudy:

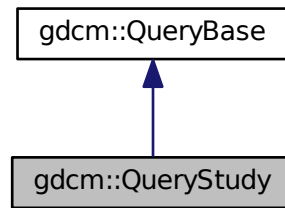

Collaboration diagram for gdcm::QueryStudy:

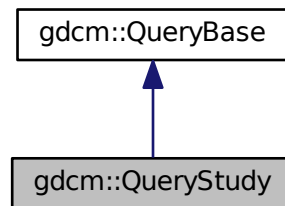

## Public Member Functions

- `std::vector< Tag > GetHierachicalSearchTags` (const `ERootType` &inRootType) const  
*Return all Unique Key for a particular Query Root type (from the same level and above).*
- `const char * GetName` () const
- `std::vector< Tag > GetOptionalTags` (const `ERootType` &inRootType) const
- `DataElement GetQueryLevel` () const
- `std::vector< Tag > GetRequiredTags` (const `ERootType` &inRootType) const
- `std::vector< Tag > GetUniqueTags` (const `ERootType` &inRootType) const

### 27.229.1 Detailed Description

QueryStudy.h contains: class to construct a study-based query for C-FIND and C-MOVE.

### 27.229.2 Member Function Documentation

**27.229.2.1** `std::vector<Tag> gdcmm::QueryStudy::GetHierachicalSearchTags ( const ERootType & inRootType ) const` `[virtual]`

Return all Unique Key for a particular Query Root type (from the same level and above).

Implements [gdcmm::QueryBase](#).

**27.229.2.2** `const char* gdcmm::QueryStudy::GetName ( ) const` `[virtual]`

Implements [gdcmm::QueryBase](#).

**27.229.2.3** `std::vector<Tag> gdcmm::QueryStudy::GetOptionalTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcmm::QueryBase](#).

**27.229.2.4** `DataElement gdcmm::QueryStudy::GetQueryLevel ( ) const` `[virtual]`

Implements [gdcmm::QueryBase](#).

**27.229.2.5** `std::vector<Tag> gdcmm::QueryStudy::GetRequiredTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcmm::QueryBase](#).

**27.229.2.6** `std::vector<Tag> gdcmm::QueryStudy::GetUniqueTags ( const ERootType & inRootType ) const` `[virtual]`

Implements [gdcmm::QueryBase](#).

The documentation for this class was generated from the following file:


- [gdcmmQueryStudy.h](#)

## 27.230 gdcmm::RAWCodec Class Reference


[RAWCodec](#) class.

```
#include <gdcmmRAWCodec.h>
```

Inheritance diagram for gdcmm::RAWCodec:

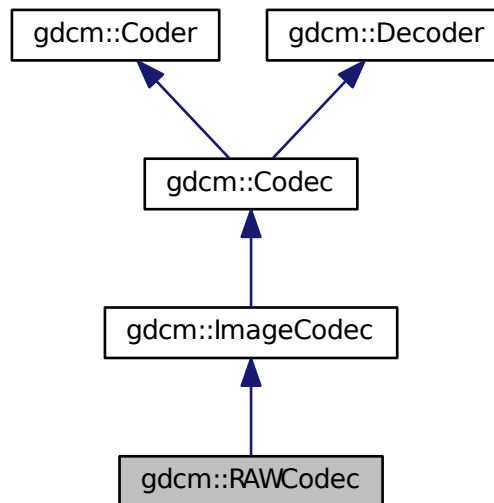

Collaboration diagram for gdcmm::RAWCodec:

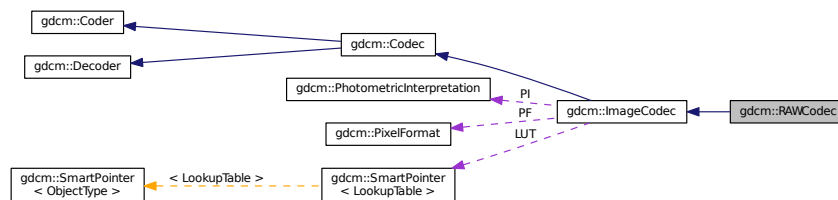

## Public Member Functions

- [RAWCodec](#) ()
- [~RAWCodec](#) ()
- bool [CanCode](#) ([TransferSyntax](#) const &ts) const  
Return whether this coder support this transfer syntax (can code it)
- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const  
Return whether this decoder support this transfer syntax (can decode it)
- virtual [ImageCodec](#) \* [Clone](#) () const
- bool [Code](#) ([DataElement](#) const &in, [DataElement](#) &out)  
Code.
- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)

*Decode.*

- bool [DecodeBytes](#) (const char \*inBytes, size\_t inBufferLength, char \*outBytes, size\_t inOutBufferLength)
- bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)

## Protected Member Functions

- bool [DecodeByStreams](#) (std::istream &is, std::ostream &os)

## Additional Inherited Members

### 27.230.1 Detailed Description

[RAWCodec](#) class.

### 27.230.2 Constructor & Destructor Documentation

27.230.2.1 `gdcm::RAWCodec::RAWCodec ( )`

27.230.2.2 `gdcm::RAWCodec::~~RAWCodec ( )`

### 27.230.3 Member Function Documentation

27.230.3.1 `bool gdcm::RAWCodec::CanCode ( TransferSyntax const & ) const` `[virtual]`

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

27.230.3.2 `bool gdcm::RAWCodec::CanDecode ( TransferSyntax const & ) const` `[virtual]`

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

27.230.3.3 `virtual ImageCodec* gdcm::RAWCodec::Clone ( ) const` `[virtual]`

Implements [gdcm::ImageCodec](#).

27.230.3.4 `bool gdcm::RAWCodec::Code ( DataElement const & in_, DataElement & out_ )` `[virtual]`

Code.

Reimplemented from [gdcm::Coder](#).

27.230.3.5 `bool gdcm::RAWCodec::Decode ( DataElement const &, DataElement & )` `[virtual]`

Decode.

Reimplemented from [gdcm::ImageCodec](#).

27.230.3.6 `bool gdcm::RAWCodec::DecodeByStreams ( std::istream & is, std::ostream & os )` `[protected]`, `[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.230.3.7 `bool gdcm::RAWCodec::DecodeBytes ( const char * inBytes, size_t inBufferLength, char * outBytes, size_t inOutBufferLength )`

Used by the ImageStreamReader— converts a read in buffer into one with the proper encodings.

27.230.3.8 `bool gdcm::RAWCodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` `[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

The documentation for this class was generated from the following file:


- [gdcmRAWCodec.h](#)

## 27.231 gdcm::Reader Class Reference


[Reader](#) ala DOM (Document [Object](#) Model)

```
#include <gdcmReader.h>
```

Inheritance diagram for `gdcm::Reader`:

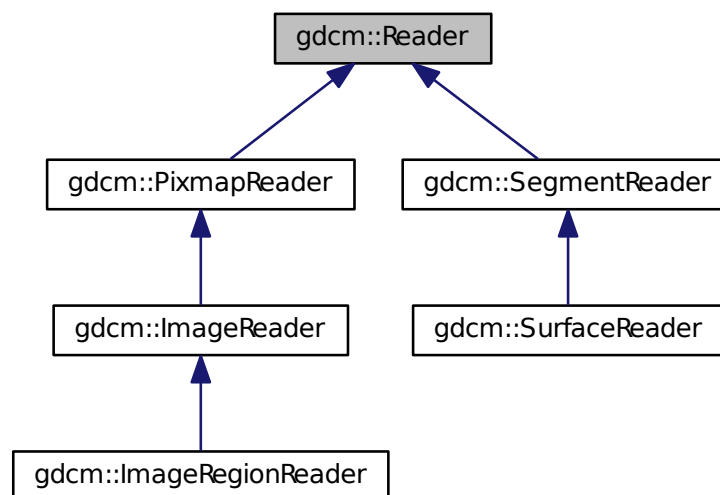

Collaboration diagram for `gdcm::Reader`:

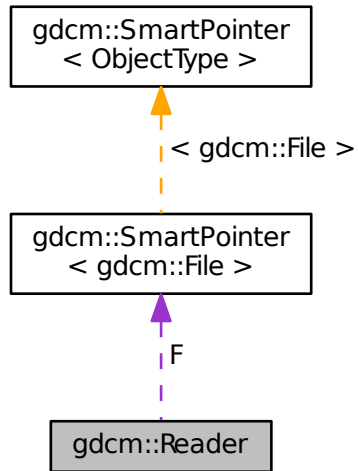

## Public Member Functions

- [Reader](#) ()
- virtual [~Reader](#) ()
- bool [CanRead](#) () const
- const [File](#) & [GetFile](#) () const  
*Set/Get File.*
- [File](#) & [GetFile](#) ()  
*Set/Get File.*
- size\_t [GetStreamCurrentPosition](#) () const
- virtual bool [Read](#) ()  
*Main function to read a file.*
- bool [ReadSelectedPrivateTags](#) (std::set< [PrivateTag](#) > const &ptags, bool readvalues=true)  
*Will only read the specified selected private tags.*
- bool [ReadSelectedTags](#) (std::set< [Tag](#) > const &tags, bool readvalues=true)  
*Will only read the specified selected tags.*
- bool [ReadUpToTag](#) (const [Tag](#) &tag, std::set< [Tag](#) > const &skiptags=std::set< [Tag](#) >())
- void [SetFile](#) ([File](#) &file)  
*Set/Get File.*
- void [SetFileName](#) (const char \*filename\_native)
- void [SetStream](#) (std::istream &input\_stream)  
*Set the open-ed stream directly.*



## Protected Member Functions

- `std::istream * GetStreamPtr () const`
- `bool ReadDataSet ()`
- `bool ReadMetaInformation ()`
- `bool ReadPreamble ()`

## Protected Attributes

- `SmartPointer< File > F`

## Friends

- `class StreamImageReader`

### 27.231.1 Detailed Description

[Reader](#) ala DOM (Document [Object](#) Model)

This class is a non-validating reader, it will only performs well- formedness check only, and to some extent catch known error (non well-formed document).

Detailed description here

A [DataSet](#) DOES NOT contains group 0x0002 (see [FileMetaInformation](#))

This is really a [DataSet](#) reader. This will not make sure the dataset conform to any [IOD](#) at all. This is a completely different step. The reasoning was that user could control the [IOD](#) there lib would handle and thus we would not be able to read a [DataSet](#) if the [IOD](#) was not found Instead we separate the reading from the validation.

#### Note

From GDCM1.x. Users will realize that one feature is missing from this DOM implementation. In GDCM 1.x user used to be able to control the size of the [Value](#) to be read. By default it was 0xffff. The main author of GDCM2 thought this was too dangerous and harmful and therefore this feature did not make it into GDCM2

#### Warning

GDCM will not produce warning for unordered (non-alphabetical order).

#### See also

[Writer](#) [FileMetaInformation](#) [DataSet](#) [File](#)

#### Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CreateFakeRTDOSE.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [Dumplmage↵HeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncrypted↵Content.cxx](#), [FixBrokenJ2K.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [Get↵SequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.↵cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [ReadUTF8QtDir.cxx](#), [rle2img.cxx](#), and [TestReader.cxx](#).

## 27.231.2 Constructor & Destructor Documentation

27.231.2.1 `gdcm::Reader::Reader ( )` `[inline]`

27.231.2.2 `virtual gdcm::Reader::~~Reader ( )` `[virtual]`

## 27.231.3 Member Function Documentation

27.231.3.1 `bool gdcm::Reader::CanRead ( )` `const`

Test whether this is a DICOM file

Warning

need to call either `SetFileName` or `SetStream` first

Examples:

[ReadUTF8QtDir.cxx](#).

27.231.3.2 `const File& gdcm::Reader::GetFile ( )` `const` `[inline]`

Set/Get [File](#).

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CompressImage.cxx](#), [CreateFakeRTDOSE.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMS↵  
MovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1Wave↵  
ToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [Get↵  
SequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [rle2img.cxx](#), and [TestReader.↵  
cxx](#).

27.231.3.3 `File& gdcm::Reader::GetFile ( )` `[inline]`

Set/Get [File](#).

27.231.3.4 `size_t gdcm::Reader::GetStreamCurrentPosition ( )` `const`

For wrapped language. return type is compatible with [System::FileSize](#) return type Use native `std::streampos` / `std↵  
::streamoff` directly from the stream from C++

27.231.3.5 `std::istream* gdcm::Reader::GetStreamPtr ( )` `const` `[inline]`, `[protected]`

27.231.3.6 `virtual bool gdcm::Reader::Read ( )` `[virtual]`

Main function to read a file.

Reimplemented in [gdcm::ImageRegionReader](#), [gdcm::PixmapReader](#), [gdcm::ImageReader](#), [gdcm::SegmentReader](#), and [gdcm::SurfaceReader](#).

#### Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CreateFakeRTDOSE.cxx](#), [csa2img.cxx](#), [DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [FixBrokenJ2K.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [rle2img.cxx](#), and [TestReader.cxx](#).

27.231.3.7 `bool gdcm::Reader::ReadDataSet ( )` [protected]

27.231.3.8 `bool gdcm::Reader::ReadMetaInformation ( )` [protected]

27.231.3.9 `bool gdcm::Reader::ReadPreamble ( )` [protected]

27.231.3.10 `bool gdcm::Reader::ReadSelectedPrivateTags ( std::set< PrivateTag > const & ptags, bool readvalues = true )`

Will only read the specified selected private tags.

27.231.3.11 `bool gdcm::Reader::ReadSelectedTags ( std::set< Tag > const & tags, bool readvalues = true )`

Will only read the specified selected tags.

27.231.3.12 `bool gdcm::Reader::ReadUpToTag ( const Tag & tag, std::set< Tag > const & skiptags = std::set< Tag > ( ) )`

Will read only up to [Tag](#)

#### Parameters

<i>tag</i>	and skipping any tag specified in
<i>skiptags</i>	

27.231.3.13 `void gdcm::Reader::SetFile ( File & file )` [inline]

Set/Get [File](#).

27.231.3.14 `void gdcm::Reader::SetFileName ( const char * filename_native )`

Set the filename to open. This will create a `std::ifstream` internally See [SetStream](#) if you are dealing with different `std::istream` object

#### Examples:

[BasicImageAnonymizer.cs](#), [ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [CheckBigEndianBug.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CompressImage.cxx](#), [ConvertToQImage.cxx](#), [CreateFakeRTDOSE.cxx](#), [csa2img.cxx](#),

[DiffFile.cxx](#), [DumpADAC.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpImageHeaderInfo.cxx](#), [DumpPhilipsECHO.cxx](#), [DuplicatePCDE.cxx](#), [ELSCINT1WaveToText.cxx](#), [ExtractEncryptedContent.cxx](#), [ExtractIconFromFile.cxx](#), [ExtractImageRegion.cs](#), [ExtractImageRegionWithLUT.cs](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [gdcmrtonplan.cxx](#), [gdcmrtpplan.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [MrProtocol.cxx](#), [NewSequence.cs](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [ReadGEMSSDO.cxx](#), [ReadMultiTimesException.cxx](#), [ReadUTF8QtDir.cxx](#), [rle2img.cxx](#), [SimplePrintPatientName.cs](#), [TestReader.cxx](#), and [threadgdcmm.cxx](#).

27.231.3.15 `void gdcmm::Reader::SetStream ( std::istream & input_stream ) [inline]`

Set the open-ed stream directly.

Examples:

[ReadUTF8QtDir.cxx](#).

## 27.231.4 Friends And Related Function Documentation

27.231.4.1 `friend class StreamImageReader [friend]`

## 27.231.5 Member Data Documentation

27.231.5.1 `SmartPointer<File> gdcmm::Reader::F [protected]`

The documentation for this class was generated from the following file:


- [gdcmmReader.h](#)

## 27.232 gdcmm::Region Class Reference

Class for manipulation region.

```
#include <gdcmmRegion.h>
```

Inheritance diagram for gdcmm::Region:

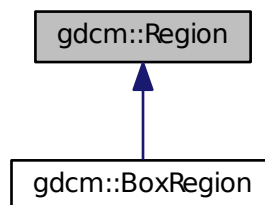

## Public Member Functions

- [Region](#) ()
- virtual [~Region](#) ()
- virtual size\_t [Area](#) () const =0  
*compute the area*
- virtual [Region](#) \* [Clone](#) () const =0
- virtual [BoxRegion](#) [ComputeBoundingBox](#) ()=0  
*Return the Axis-Aligned minimum bounding box for all regions.*
- virtual bool [Empty](#) () const =0  
*return whether this domain is empty:*
- virtual bool [IsValid](#) () const =0  
*return whether this is valid domain*
- virtual void [Print](#) (std::ostream &os=std::cout) const  
*Print.*

### 27.232.1 Detailed Description

Class for manipulation region.

### 27.232.2 Constructor & Destructor Documentation

27.232.2.1 `gdcm::Region::Region ( )`

27.232.2.2 `virtual gdcm::Region::~~Region ( )` [virtual]

### 27.232.3 Member Function Documentation

27.232.3.1 `virtual size_t gdcm::Region::Area ( ) const` [pure virtual]

compute the area

Implemented in [gdcm::BoxRegion](#).

27.232.3.2 `virtual Region* gdcm::Region::Clone ( ) const` [pure virtual]

Implemented in [gdcm::BoxRegion](#).

27.232.3.3 `virtual BoxRegion gdcm::Region::ComputeBoundingBox ( )` [pure virtual]

Return the Axis-Aligned minimum bounding box for all regions.

Implemented in [gdcm::BoxRegion](#).

27.232.3.4 `virtual bool gdcm::Region::Empty ( ) const` [pure virtual]

return whether this domain is empty:

Implemented in [gdcm::BoxRegion](#).

27.232.3.5 `virtual bool gdcmm::Region::IsValid ( ) const [pure virtual]`

return whether this is valid domain

Implemented in [gdcmm::BoxRegion](#).

27.232.3.6 `virtual void gdcmm::Region::Print ( std::ostream & os = std::cout ) const [virtual]`

Print.

Reimplemented in [gdcmm::BoxRegion](#).

Referenced by `gdcmm::operator<<()`.

The documentation for this class was generated from the following file:

- [gdcmmRegion.h](#)

## 27.233 gdcmm::Rescaler Class Reference

**Rescale class** This class is meant to apply the linear transform of Stored Pixel [Value](#) to Real World [Value](#). This is mostly found in CT or PET dataset, where the value are stored using one type, but need to be converted to another scale using a linear transform. There are basically two cases: In CT: the linear transform is generally integer based. E.g. the Stored Pixel [Type](#) is unsigned short 12bits, but to get Hounsfield unit, one need to apply the linear transform:

$$RWV = 1. * SV - 1024$$

So the best scalar to store the Real World [Value](#) will be 16 bits signed type.

```
#include <gdcmmRescaler.h>
```

### Public Member Functions

- [Rescaler](#) ()
- [~Rescaler](#) ()
- [PixelFormat::ScalarType ComputeInterceptSlopePixelFormat](#) ()
- [PixelFormat ComputePixelFormatFromMinMax](#) ()
- double [GetIntercept](#) () const
- double [GetSlope](#) () const
- bool [InverseRescale](#) (char \*out, const char \*in, size\_t n)  
*Inverse transform.*
- bool [Rescale](#) (char \*out, const char \*in, size\_t n)  
*Direct transform.*
- void [SetIntercept](#) (double i)  
*Set Intercept: used for both direct&inverse transformation.*
- void [SetMinMaxForPixelFormat](#) (double min, double max)
- void [SetPixelFormat](#) ([PixelFormat](#) const &pf)  
*Set Pixel Format of input data.*
- void [SetSlope](#) (double s)  
*Set Slope: user for both direct&inverse transformation.*
- void [SetTargetPixelFormat](#) ([PixelFormat](#) const &targetst)
- void [SetUseTargetPixelFormat](#) (bool b)  
*Override default behavior of Rescale.*

## Protected Member Functions

- template<typename TIn >  
void [InverseRescaleFunctionIntoBestFit](#) (char \*out, const TIn \*in, size\_t n)
- template<typename TIn >  
void [RescaleFunctionIntoBestFit](#) (char \*out, const TIn \*in, size\_t n)

### 27.233.1 Detailed Description

Rescale class This class is meant to apply the linear transform of Stored Pixel [Value](#) to Real World [Value](#). This is mostly found in CT or PET dataset, where the value are stored using one type, but need to be converted to another scale using a linear transform. There are basically two cases: In CT: the linear transform is generally integer based. E.g. the Stored Pixel [Type](#) is unsigned short 12bits, but to get Hounsfield unit, one need to apply the linear transform:

$$RWV = 1. * SV - 1024$$

So the best scalar to store the Real World [Value](#) will be 16 bits signed type.

In PET: the linear transform is generally floating point based. Since the dynamic range can be quite high, the Rescale Slope / Rescale Intercept can be changing throughout the [Series](#). So it is important to read all linear transform and deduce the best Pixel [Type](#) only at the end (when all the images to be read have been parsed).

#### Warning

Internally any time a floating point value is found either in the Rescale Slope or the Rescale Intercept it is assumed that the best matching output pixel type is FLOAT64 (in previous implementation it was FLOAT32). Because [V↔R:DS](#) is closer to a 64bits floating point type FLOAT64 is thus a best matching pixel type for the floating point transformation.

Example: Let say input is FLOAT64, and we want UINT16 as ouput, we would do:

```
Rescaler ir;
ir.SetIntercept( 0 );
ir.SetSlope( 5.6789 );
ir.SetPixelFormat( FLOAT64 );
ir.SetMinMaxForPixelType( ((PixelFormat)UINT16).GetMin(), ((PixelFormat)UINT16).GetMax() );
ir.InverseRescale(output,input,numberofbytes );
```

#### Note

handle floating point transformation back and forth to integer properly (no loss)

#### See also

[Unpacker12Bits](#)

### 27.233.2 Constructor & Destructor Documentation

27.233.2.1 `gdcm::Rescaler::Rescaler ( )` `[inline]`

27.233.2.2 `gdcm::Rescaler::~~Rescaler ( )` `[inline]`

### 27.233.3 Member Function Documentation

27.233.3.1 `PixelFormat::ScalarType gdcm::Rescaler::ComputeInterceptSlopePixelFormat ( )`

Compute the Pixel Format of the output data Used for direct transformation

### 27.233.3.2 PixelFormat gdcm::Rescaler::ComputePixelTypeFromMinMax ( )

Compute the Pixel Format of the output data Used for inverse transformation

### 27.233.3.3 double gdcm::Rescaler::GetIntercept ( ) const [inline]

### 27.233.3.4 double gdcm::Rescaler::GetSlope ( ) const [inline]

### 27.233.3.5 bool gdcm::Rescaler::InverseRescale ( char \* out, const char \* in, size\_t n )

Inverse transform.

### 27.233.3.6 template<typename TIn > void gdcm::Rescaler::InverseRescaleFunctionIntoBestFit ( char \* out, const TIn \* in, size\_t n ) [protected]

### 27.233.3.7 bool gdcm::Rescaler::Rescale ( char \* out, const char \* in, size\_t n )

Direct transform.

### 27.233.3.8 template<typename TIn > void gdcm::Rescaler::RescaleFunctionIntoBestFit ( char \* out, const TIn \* in, size\_t n ) [protected]

### 27.233.3.9 void gdcm::Rescaler::SetIntercept ( double i ) [inline]

Set Intercept: used for both direct&inverse transformation.

### 27.233.3.10 void gdcm::Rescaler::SetMinMaxForPixelType ( double min, double max ) [inline]

Set target interval for output data. A best match will be computed (if possible) Used for inverse transformation

### 27.233.3.11 void gdcm::Rescaler::SetPixelFormat ( PixelFormat const & pf ) [inline]

Set Pixel Format of input data.

### 27.233.3.12 void gdcm::Rescaler::SetSlope ( double s ) [inline]

Set Slope: user for both direct&inverse transformation.

### 27.233.3.13 void gdcm::Rescaler::SetTargetPixelType ( PixelFormat const & targetst )

By default (when UseTargetPixelType is false), a best matching Target Pixel [Type](#) is computed. However user can override this auto selection by switching UseTargetPixelType:true and also specifying the specifix Target Pixel [Type](#)

### 27.233.3.14 void gdcm::Rescaler::SetUseTargetPixelType ( bool b )

Override default behavior of Rescale.

The documentation for this class was generated from the following file:




- [gdcmRescaler.h](#)

## 27.234 gdcm::RLECodec Class Reference


Class to do RLE.

```
#include <gdcmRLECodec.h>
```

Inheritance diagram for gdcm::RLECodec:

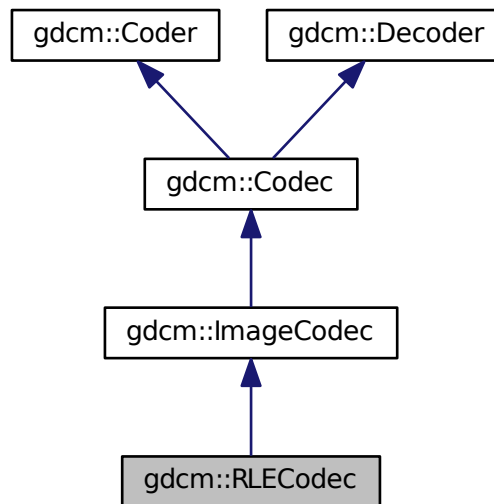

Collaboration diagram for gdcm::RLECodec:

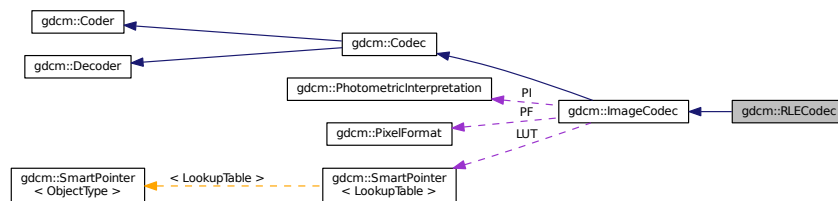

### Public Member Functions

- [RLECodec](#) ()
- [~RLECodec](#) ()
- [bool CanCode](#) ([TransferSyntax](#) const &ts) const

*Return whether this coder support this transfer syntax (can code it)*

- bool [CanDecode](#) ([TransferSyntax](#) const &ts) const

*Return whether this decoder support this transfer syntax (can decode it)*

- virtual [ImageCodec](#) \* [Clone](#) () const
- bool [Code](#) ([DataElement](#) const &in, [DataElement](#) &out)

*Code.*

- bool [Decode](#) ([DataElement](#) const &is, [DataElement](#) &os)

*Decode.*

- unsigned long [GetBufferLength](#) () const
- bool [GetHeaderInfo](#) (std::istream &is, [TransferSyntax](#) &ts)
- void [SetBufferLength](#) (unsigned long l)
- void [SetLength](#) (unsigned long l)

## Protected Member Functions

- bool [AppendFrameEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [AppendRowEncode](#) (std::ostream &out, const char \*data, size\_t datalen)
- bool [DecodeByStreams](#) (std::istream &is, std::ostream &os)
- bool [DecodeExtent](#) (char \*buffer, unsigned int XMin, unsigned int XMax, unsigned int YMin, unsigned int YMax, unsigned int ZMin, unsigned int ZMax, std::istream &is)
- bool [IsFrameEncoder](#) ()
- bool [IsRowEncoder](#) ()
- bool [StartEncode](#) (std::ostream &)
- bool [StopEncode](#) (std::ostream &)

## Friends

- class [ImageRegionReader](#)

## Additional Inherited Members

### 27.234.1 Detailed Description

Class to do RLE.

#### Note

ANSI X3.9 A.4.2 RLE Compression Annex G defines a RLE Compression Transfer Syntax. This transfer Syntax is identified by the UID value "1.2.840.10008.1.2.5". If the object allows multi-frame images in the pixel data field, then each frame shall be encoded separately. Each frame shall be encoded in one and only one [Fragment](#) (see PS 3.5.8.2).

### 27.234.2 Constructor & Destructor Documentation

27.234.2.1 [gdcm::RLECodec::RLECodec](#) ( )

27.234.2.2 [gdcm::RLECodec::~~RLECodec](#) ( )

### 27.234.3 Member Function Documentation

**27.234.3.1** `bool gdcm::RLECodec::AppendFrameEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.2** `bool gdcm::RLECodec::AppendRowEncode ( std::ostream & out, const char * data, size_t datalen )`  
[protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.3** `bool gdcm::RLECodec::CanCode ( TransferSyntax const & ) const` [virtual]

Return whether this coder support this transfer syntax (can code it)

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.4** `bool gdcm::RLECodec::CanDecode ( TransferSyntax const & ) const` [virtual]

Return whether this decoder support this transfer syntax (can decode it)

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.5** `virtual ImageCodec* gdcm::RLECodec::Clone ( ) const` [virtual]

Implements [gdcm::ImageCodec](#).

**27.234.3.6** `bool gdcm::RLECodec::Code ( DataElement const & in_, DataElement & out_ )` [virtual]

Code.

Reimplemented from [gdcm::Coder](#).

**27.234.3.7** `bool gdcm::RLECodec::Decode ( DataElement const &, DataElement & )` [virtual]

Decode.

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.8** `bool gdcm::RLECodec::DecodeByStreams ( std::istream & is, std::ostream & os )` [protected], [virtual]

Reimplemented from [gdcm::ImageCodec](#).

**27.234.3.9** `bool gdcm::RLECodec::DecodeExtent ( char * buffer, unsigned int XMin, unsigned int XMax, unsigned int YMin, unsigned int YMax, unsigned int ZMin, unsigned int ZMax, std::istream & is )` [protected]

**27.234.3.10** `unsigned long gdcm::RLECodec::GetBufferLength ( ) const` [inline]

**27.234.3.11** `bool gdcm::RLECodec::GetHeaderInfo ( std::istream & is, TransferSyntax & ts )` [virtual]

Reimplemented from [gdcm::ImageCodec](#).

27.234.3.12 `bool gdcm::RLECodec::IsFrameEncoder ( ) [protected],[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.234.3.13 `bool gdcm::RLECodec::IsRowEncoder ( ) [protected],[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.234.3.14 `void gdcm::RLECodec::SetBufferLength ( unsigned long l ) [inline]`

27.234.3.15 `void gdcm::RLECodec::SetLength ( unsigned long l ) [inline]`

27.234.3.16 `bool gdcm::RLECodec::StartEncode ( std::ostream & ) [protected],[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

27.234.3.17 `bool gdcm::RLECodec::StopEncode ( std::ostream & ) [protected],[virtual]`

Reimplemented from [gdcm::ImageCodec](#).

## 27.234.4 Friends And Related Function Documentation

27.234.4.1 `friend class ImageRegionReader [friend]`

The documentation for this class was generated from the following file:

- [gdcmRLECodec.h](#)

## 27.235 gdcm::network::RoleSelectionSub Class Reference

[RoleSelectionSub](#) PS 3.7 [Table D.3-9](#) SCP/SCU ROLE SELECTION SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

```
#include <gdcmRoleSelectionSub.h>
```

### Public Member Functions

- [RoleSelectionSub](#) ( )
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetTuple](#) (const char \*uid, uint8\_t scurole, uint8\_t scprole)
- size\_t [Size](#) ( ) const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.235.1 Detailed Description

[RoleSelectionSub](#) PS 3.7 [Table D.3-9](#) SCP/SCU ROLE SELECTION SUB-ITEM FIELDS (A-ASSOCIATE-RQ)

## 27.235.2 Constructor & Destructor Documentation

27.235.2.1 `gdcm::network::RoleSelectionSub::RoleSelectionSub ( )`

## 27.235.3 Member Function Documentation

27.235.3.1 `void gdcm::network::RoleSelectionSub::Print ( std::ostream & os ) const`

27.235.3.2 `std::istream& gdcm::network::RoleSelectionSub::Read ( std::istream & is )`

27.235.3.3 `void gdcm::network::RoleSelectionSub::SetTuple ( const char * uid, uint8_t scurole, uint8_t scprole )`

27.235.3.4 `size_t gdcm::network::RoleSelectionSub::Size ( ) const`

27.235.3.5 `const std::ostream& gdcm::network::RoleSelectionSub::Write ( std::ostream & os ) const`


The documentation for this class was generated from the following file:

- [gdcmRoleSelectionSub.h](#)

## 27.236 gdcm::SerieHelper::Rule Struct Reference

```
#include <gdcmSerieHelper.h>
```

Collaboration diagram for `gdcm::SerieHelper::Rule`:

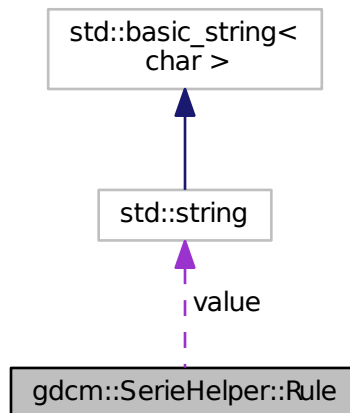

## Public Attributes

- `uint16_t elem`
- `uint16_t group`

- int [op](#)
- std::string [value](#)

### 27.236.1 Member Data Documentation

27.236.1.1 uint16\_t [gdcm::SerieHelper::Rule::elem](#)

27.236.1.2 uint16\_t [gdcm::SerieHelper::Rule::group](#)

27.236.1.3 int [gdcm::SerieHelper::Rule::op](#)

27.236.1.4 std::string [gdcm::SerieHelper::Rule::value](#)

The documentation for this struct was generated from the following file:


- [gdcmSerieHelper.h](#)

## 27.237 gdcm::Scanner Class Reference


[Scanner](#) This filter is meant for quickly browsing a [FileSet](#) (a set of files on disk). Special consideration are taken so as to read the minimum amount of information in each file in order to retrieve the user specified set of DICOM [Attribute](#).

```
#include <gdcmScanner.h>
```

Inheritance diagram for [gdcm::Scanner](#):

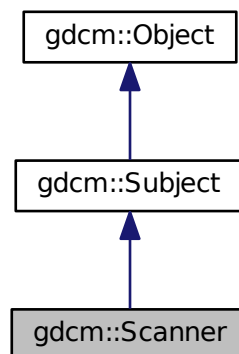

Collaboration diagram for gdcm::Scanner:

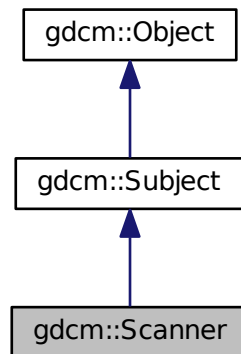

## Classes

- struct [Itstr](#)

## Public Types

- typedef MappingType::const\_iterator [ConstIterator](#)
- typedef std::map< const char \*, [TagToValue](#), [Itstr](#) > [MappingType](#)
- typedef std::map< [Tag](#), const char \* > [TagToValue](#)
- typedef TagToValue::value\_type [TagToValueValueType](#)
- typedef std::set< std::string > [ValuesType](#)

## Public Member Functions

- [Scanner](#) ()
- [~Scanner](#) ()
- void [AddPrivateTag](#) ([PrivateTag](#) const &t)
- void [AddSkipTag](#) ([Tag](#) const &t)  
*Add a tag that will need to be skipped. Those are root level skip tags.*
- void [AddTag](#) ([Tag](#) const &t)  
*Add a tag that will need to be read. Those are root level skip tags.*
- [ConstIterator](#) [Begin](#) () const
- void [ClearSkipTags](#) ()
- void [ClearTags](#) ()
- [ConstIterator](#) [End](#) () const
- [Directory::FilenameType](#) [GetAllFileNamesFromTagToValue](#) ([Tag](#) const &t, const char \*valueref) const
- const char \* [GetFilenameFromTagToValue](#) ([Tag](#) const &t, const char \*valueref) const
- [Directory::FilenameType](#) const & [GetFileNames](#) () const

- [Directory::FilenamesType GetKeys](#) () const
- [TagToValue](#) const & [GetMapping](#) (const char \*filename) const  
*Get the std::map mapping filenames to value for file 'filename'.*
- [TagToValue](#) const & [GetMappingFromTagToValue](#) ([Tag](#) const &t, const char \*value) const  
*See [GetFilenameFromTagToValue\(\)](#). This is simply [GetFilenameFromTagToValue](#) followed.*
- [MappingType](#) const & [GetMappings](#) () const  
*Mappings are the mapping from a particular tag to the map, mapping filename to value:*
- [Directory::FilenamesType GetOrderedValues](#) ([Tag](#) const &t) const
- const char \* [GetValue](#) (const char \*filename, [Tag](#) const &t) const
- [ValuesType](#) const & [GetValues](#) () const  
*Get all the values found (in lexicographic order)*
- [ValuesType GetValues](#) ([Tag](#) const &t) const  
*Get all the values found (in lexicographic order) associated with [Tag](#) 't'.*
- bool [IsKey](#) (const char \*filename) const
- void [Print](#) (std::ostream &os) const  
*Print result.*
- bool [Scan](#) ([Directory::FilenamesType](#) const &filenames)  
*Start the scan !*

## Static Public Member Functions

- static [SmartPointer](#)< [Scanner](#) > [New](#) ()  
*for wrapped language: instantiate a reference counted object*

## Protected Member Functions

- void [ProcessPublicTag](#) ([StringFilter](#) &sf, const char \*filename)

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Scanner](#) &s)

## 27.237.1 Detailed Description

[Scanner](#) This filter is meant for quickly browsing a [FileSet](#) (a set of files on disk). Special consideration are taken so as to read the minimum amount of information in each file in order to retrieve the user specified set of DICOM [Attribute](#).

This filter is dealing with both VRASCII and VRBINARY element, thanks to the help of [StringFilter](#)

### Warning

IMPORTANT In case of file where tags are not ordered (illegal as per DICOM specification), the output will be missing information



## Note

implementation details. All values are stored in a `std::set` of `std::string`. Then the address of the `cstring` underlying the `std::string` is used in the `std::map`.

This class implement the Subject/Observer pattern trigger the following events:

- [ProgressEvent](#)
- [StartEvent](#)
- [EndEvent](#)

## Examples:

[DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [SimpleScanner.cxx](#), [SortImage.cxx](#), and [VolumeSorter.cxx](#).

## 27.237.2 Member Typedef Documentation

27.237.2.1 `typedef MappingType::const_iterator gdcm::Scanner::ConstIterator`

27.237.2.2 `typedef std::map<const char *, TagToValue, Itstr> gdcm::Scanner::MappingType`

27.237.2.3 `typedef std::map<Tag, const char*> gdcm::Scanner::TagToValue`

struct to map a filename to a value Implementation note: all `std::map` in this class will be using `const char *` and not `std::string` since we are pointing to existing `std::string` (hold in a `std::vector`) this avoid an extra copy of the byte array. `Tag` are used as `Tag` class since `sizeof(tag) <= sizeof(pointer)`

27.237.2.4 `typedef TagToValue::value_type gdcm::Scanner::TagToValueValueType`

27.237.2.5 `typedef std::set< std::string > gdcm::Scanner::ValuesType`

## 27.237.3 Constructor &amp; Destructor Documentation

27.237.3.1 `gdcm::Scanner::Scanner ( ) [inline]`

27.237.3.2 `gdcm::Scanner::~~Scanner ( )`

## 27.237.4 Member Function Documentation

27.237.4.1 `void gdcm::Scanner::AddPrivateTag ( PrivateTag const & t )`

27.237.4.2 `void gdcm::Scanner::AddSkipTag ( Tag const & t )`

Add a tag that will need to be skipped. Those are root level skip tags.

27.237.4.3 `void gdcm::Scanner::AddTag ( Tag const & t )`

Add a tag that will need to be read. Those are root level skip tags.

## Examples:

[DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [SimpleScanner.cxx](#), [SortImage.cxx](#), and [VolumeSorter.cxx](#).

27.237.4.4 **ConstIterator** `gdcmm::Scanner::Begin ( ) const` `[inline]`

27.237.4.5 `void gdcmm::Scanner::ClearSkipTags ( )`

27.237.4.6 `void gdcmm::Scanner::ClearTags ( )`

27.237.4.7 **ConstIterator** `gdcmm::Scanner::End ( ) const` `[inline]`

27.237.4.8 **Directory::FilenameType** `gdcmm::Scanner::GetAllFilenamesFromTagToValue ( Tag const & t, const char * valuref ) const`

Will loop over all files and return a vector of `std::strings` of filenames where value match the reference value 'valuref'

27.237.4.9 `const char* gdcmm::Scanner::GetFilenameFromTagToValue ( Tag const & t, const char * valuref ) const`

Will loop over all files and return the first file where value match the reference value 'valuref'

27.237.4.10 **Directory::FilenameType** `const& gdcmm::Scanner::GetFilenames ( ) const` `[inline]`

27.237.4.11 **Directory::FilenameType** `gdcmm::Scanner::GetKeys ( ) const`

Return the list of filename that are key in the internal map, which means those filename were properly parsed

Examples:

[VolumeSorter.cxx](#).

27.237.4.12 **TagToValue** `const& gdcmm::Scanner::GetMapping ( const char * filename ) const`

Get the `std::map` mapping filenames to value for file 'filename'.

Examples:

[DumpToSQLITE3.cxx](#), and [SimpleScanner.cxx](#).

27.237.4.13 **TagToValue** `const& gdcmm::Scanner::GetMappingFromTagToValue ( Tag const & t, const char * value ) const`

See [GetFilenameFromTagToValue\(\)](#). This is simply `GetFilenameFromTagToValue` followed.

27.237.4.14 **MappingType** `const& gdcmm::Scanner::GetMappings ( ) const` `[inline]`

Mappings are the mapping from a particular tag to the map, mapping filename to value:

27.237.4.15 **Directory::FilenameType** `gdcmm::Scanner::GetOrderedValues ( Tag const & t ) const`

Get all the values found (in a vector) associated with [Tag](#) 't' This function is identical to `GetValues`, but is accessible from the wrapped layer (python, C#, java)

**27.237.4.16** `const char* gdcm::Scanner::GetValue ( const char * filename, Tag const & t ) const`

Retrieve the value found for tag: t associated with file: filename This is meant for a single short call. If multiple calls (multiple tags) should be done, prefer the GetMapping function, and then reuse the TagToValue hash table.

Warning

[Tag](#) 't' should have been added via [AddTag\(\)](#) prior to the [Scan\(\)](#) call !

**27.237.4.17** `ValueType const& gdcm::Scanner::GetValues ( ) const [inline]`

Get all the values found (in lexicographic order)

Examples:

[SortImage.cxx](#), and [VolumeSorter.cxx](#).

**27.237.4.18** `ValueType gdcm::Scanner::GetValues ( Tag const & t ) const`

Get all the values found (in lexicographic order) associated with [Tag](#) 't'.

**27.237.4.19** `bool gdcm::Scanner::IsKey ( const char * filename ) const`

Check if filename is a key in the Mapping table. returns true only if file can be found, which means the file was indeed a DICOM file that could be processed

Examples:

[DumpToSQLITE3.cxx](#), and [SimpleScanner.cxx](#).

**27.237.4.20** `static SmartPointer<Scanner> gdcm::Scanner::New ( ) [inline],[static]`

for wrapped language: instantiate a reference counted object

**27.237.4.21** `void gdcm::Scanner::Print ( std::ostream & os ) const [virtual]`

Print result.

Reimplemented from [gdcm::Object](#).

Referenced by `gdcm::operator<<()`.

**27.237.4.22** `void gdcm::Scanner::ProcessPublicTag ( StringFilter & sf, const char * filename ) [protected]`

**27.237.4.23** `bool gdcm::Scanner::Scan ( Directory::FileNamesType const & filenames )`

Start the scan !

Examples:

[DiscriminateVolume.cxx](#), [DumpToSQLITE3.cxx](#), [SimpleScanner.cxx](#), [SortImage.cxx](#), and [VolumeSorter.cxx](#).

## 27.237.5 Friends And Related Function Documentation

27.237.5.1 `std::ostream& operator<< ( std::ostream & _os, const Scanner & s )` [*friend*]

The documentation for this class was generated from the following file:


- [gdcmscanner.h](#)

## 27.238 gdcmscanner::Segment Class Reference


This class defines a segment. It mainly contains attributes of group 0x0062. In addition, it can be associated with surface.

```
#include <gdcmscanner.h>
```

Inheritance diagram for `gdcmscanner::Segment`:

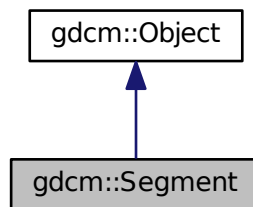

Collaboration diagram for `gdcmscanner::Segment`:

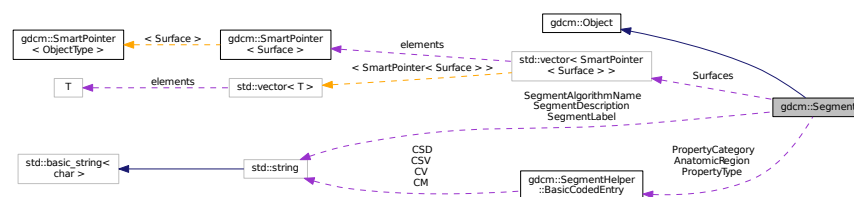

## Public Types

- enum `ALGOType` {  
`MANUAL` = 0,  
`AUTOMATIC`,  
`ALGOType_END` }
- typedef `std::vector< SmartPointer< Surface > >` `SurfaceVector`

## Public Member Functions

- [Segment](#) ()
- virtual [~Segment](#) ()
- void [AddSurface](#) ([SmartPointer](#)< [Surface](#) > surface)
- [SegmentHelper::BasicCodedEntry](#) const & [GetAnatomicRegion](#) () const
- [SegmentHelper::BasicCodedEntry](#) & [GetAnatomicRegion](#) ()
- [SegmentHelper::BasicCodedEntry](#) const & [GetPropertyCategory](#) () const
- [SegmentHelper::BasicCodedEntry](#) & [GetPropertyCategory](#) ()
- [SegmentHelper::BasicCodedEntry](#) const & [GetPropertyType](#) () const
- [SegmentHelper::BasicCodedEntry](#) & [GetPropertyType](#) ()
- const char \* [GetSegmentAlgorithmName](#) () const
- [ALGOType](#) [GetSegmentAlgorithmType](#) () const
- const char \* [GetSegmentDescription](#) () const
- const char \* [GetSegmentLabel](#) () const
- unsigned short [GetSegmentNumber](#) () const
- [SmartPointer](#)< [Surface](#) > [GetSurface](#) (const unsigned int idx=0) const
- unsigned long [GetSurfaceCount](#) ()
- [SurfaceVector](#) const & [GetSurfaces](#) () const
- [SurfaceVector](#) & [GetSurfaces](#) ()
- void [SetAnatomicRegion](#) ([SegmentHelper::BasicCodedEntry](#) const &BSE)
- void [SetPropertyCategory](#) ([SegmentHelper::BasicCodedEntry](#) const &BSE)
- void [SetPropertyType](#) ([SegmentHelper::BasicCodedEntry](#) const &BSE)
- void [SetSegmentAlgorithmName](#) (const char \*name)
- void [SetSegmentAlgorithmType](#) ([ALGOType](#) type)
- void [SetSegmentAlgorithmType](#) (const char \*typeStr)
- void [SetSegmentDescription](#) (const char \*description)
- void [SetSegmentLabel](#) (const char \*label)
- void [SetSegmentNumber](#) (const unsigned short num)
- void [SetSurfaceCount](#) (const unsigned long nb)

## Static Public Member Functions

- static [ALGOType](#) [GetALGOType](#) (const char \*type)
- static const char \* [GetALGOTypeString](#) ([ALGOType](#) type)

## Protected Attributes

- [SegmentHelper::BasicCodedEntry](#) [AnatomicRegion](#)
- [SegmentHelper::BasicCodedEntry](#) [PropertyCategory](#)
- [SegmentHelper::BasicCodedEntry](#) [PropertyType](#)
- std::string [SegmentAlgorithmName](#)
- [ALGOType](#) [SegmentAlgorithmType](#)
- std::string [SegmentDescription](#)
- std::string [SegmentLabel](#)
- unsigned short [SegmentNumber](#)
- unsigned long [SurfaceCount](#)
- [SurfaceVector](#) [Surfaces](#)

## Additional Inherited Members

### 27.238.1 Detailed Description

This class defines a segment. It mainly contains attributes of group 0x0062. In addition, it can be associated with surface.

See also

PS 3.3 C.8.20.2 and C.8.23

### 27.238.2 Member Typedef Documentation

27.238.2.1 `typedef std::vector< SmartPointer< Surface > > gdcm::Segment::SurfaceVector`

### 27.238.3 Member Enumeration Documentation

27.238.3.1 `enum gdcm::Segment::ALGOType`

Enumerator

***MANUAL***

***AUTOMATIC***

***ALGOType\_END***

### 27.238.4 Constructor & Destructor Documentation

27.238.4.1 `gdcm::Segment::Segment ( )`

27.238.4.2 `virtual gdcm::Segment::~~Segment ( ) [virtual]`

### 27.238.5 Member Function Documentation

27.238.5.1 `void gdcm::Segment::AddSurface ( SmartPointer< Surface > surface )`

27.238.5.2 `static ALGOType gdcm::Segment::GetALGOType ( const char * type ) [static]`

27.238.5.3 `static const char* gdcm::Segment::GetALGOTypeString ( ALGOType type ) [static]`

27.238.5.4 `SegmentHelper::BasicCodedEntry const& gdcm::Segment::GetAnatomicRegion ( ) const`

27.238.5.5 `SegmentHelper::BasicCodedEntry& gdcm::Segment::GetAnatomicRegion ( )`

27.238.5.6 `SegmentHelper::BasicCodedEntry const& gdcm::Segment::GetPropertyCategory ( ) const`

27.238.5.7 `SegmentHelper::BasicCodedEntry& gdcm::Segment::GetPropertyCategory ( )`

27.238.5.8 `SegmentHelper::BasicCodedEntry const& gdcm::Segment::GetPropertyType ( ) const`

27.238.5.9 `SegmentHelper::BasicCodedEntry& gdcm::Segment::GetPropertyType ( )`

- 27.238.5.10 `const char* gdcm::Segment::GetSegmentAlgorithmName ( ) const`
- 27.238.5.11 `ALGOType gdcm::Segment::GetSegmentAlgorithmType ( ) const`
- 27.238.5.12 `const char* gdcm::Segment::GetSegmentDescription ( ) const`
- 27.238.5.13 `const char* gdcm::Segment::GetSegmentLabel ( ) const`
- 27.238.5.14 `unsigned short gdcm::Segment::GetSegmentNumber ( ) const`
- 27.238.5.15 `SmartPointer< Surface > gdcm::Segment::GetSurface ( const unsigned int idx = 0 ) const`
- 27.238.5.16 `unsigned long gdcm::Segment::GetSurfaceCount ( )`
- 27.238.5.17 `SurfaceVector const& gdcm::Segment::GetSurfaces ( ) const`
- 27.238.5.18 `SurfaceVector& gdcm::Segment::GetSurfaces ( )`
- 27.238.5.19 `void gdcm::Segment::SetAnatomicRegion ( SegmentHelper::BasicCodedEntry const & BSE )`
- 27.238.5.20 `void gdcm::Segment::SetPropertyCategory ( SegmentHelper::BasicCodedEntry const & BSE )`
- 27.238.5.21 `void gdcm::Segment::SetPropertyType ( SegmentHelper::BasicCodedEntry const & BSE )`
- 27.238.5.22 `void gdcm::Segment::SetSegmentAlgorithmName ( const char * name )`
- 27.238.5.23 `void gdcm::Segment::SetSegmentAlgorithmType ( ALGOType type )`
- 27.238.5.24 `void gdcm::Segment::SetSegmentAlgorithmType ( const char * typeStr )`
- 27.238.5.25 `void gdcm::Segment::SetSegmentDescription ( const char * description )`
- 27.238.5.26 `void gdcm::Segment::SetSegmentLabel ( const char * label )`
- 27.238.5.27 `void gdcm::Segment::SetSegmentNumber ( const unsigned short num )`
- 27.238.5.28 `void gdcm::Segment::SetSurfaceCount ( const unsigned long nb )`

## 27.238.6 Member Data Documentation

- 27.238.6.1 `SegmentHelper::BasicCodedEntry gdcm::Segment::AnatomicRegion` [protected]
- 27.238.6.2 `SegmentHelper::BasicCodedEntry gdcm::Segment::PropertyCategory` [protected]
- 27.238.6.3 `SegmentHelper::BasicCodedEntry gdcm::Segment::PropertyType` [protected]
- 27.238.6.4 `std::string gdcm::Segment::SegmentAlgorithmName` [protected]
- 27.238.6.5 `ALGOType gdcm::Segment::SegmentAlgorithmType` [protected]
- 27.238.6.6 `std::string gdcm::Segment::SegmentDescription` [protected]

27.238.6.7 `std::string gdcM::Segment::SegmentLabel` [protected]

27.238.6.8 `unsigned short gdcM::Segment::SegmentNumber` [protected]

27.238.6.9 `unsigned long gdcM::Segment::SurfaceCount` [protected]

27.238.6.10 `SurfaceVector gdcM::Segment::Surfaces` [protected]

The documentation for this class was generated from the following file:


- [gdcMSegment.h](#)

## 27.239 gdcM::SegmentedPaletteColorLookupTable Class Reference

[SegmentedPaletteColorLookupTable](#) class.


```
#include <gdcMSegmentedPaletteColorLookupTable.h>
```

Inheritance diagram for gdcM::SegmentedPaletteColorLookupTable:

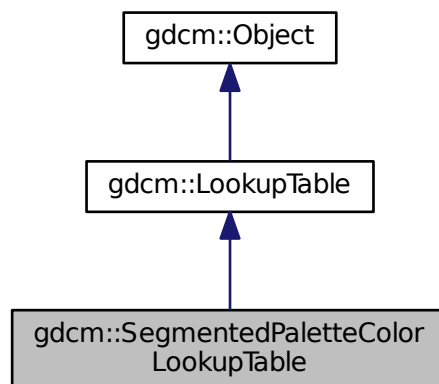



Collaboration diagram for gdcm::SegmentedPaletteColorLookupTable:

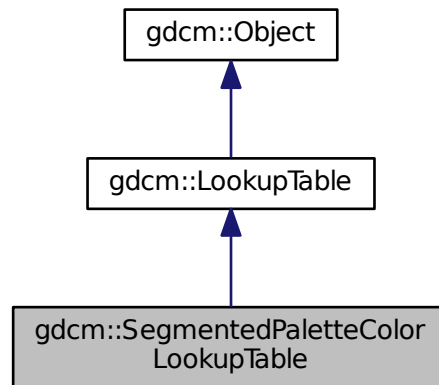

## Public Member Functions

- [SegmentedPaletteColorLookupTable](#) ()
- [~SegmentedPaletteColorLookupTable](#) ()
- void [Print](#) (std::ostream &) const
- void [SetLUT](#) ([LookupTableType](#) type, const unsigned char \*array, unsigned int length)

*Initialize a [SegmentedPaletteColorLookupTable](#).*

## Additional Inherited Members

### 27.239.1 Detailed Description

[SegmentedPaletteColorLookupTable](#) class.

### 27.239.2 Constructor & Destructor Documentation

27.239.2.1 `gdcm::SegmentedPaletteColorLookupTable::SegmentedPaletteColorLookupTable ( )`

27.239.2.2 `gdcm::SegmentedPaletteColorLookupTable::~~SegmentedPaletteColorLookupTable ( )`

### 27.239.3 Member Function Documentation

27.239.3.1 `void gdcm::SegmentedPaletteColorLookupTable::Print ( std::ostream & ) const` `[inline]`, `[virtual]`

Reimplemented from [gdcm::LookupTable](#).

27.239.3.2 void `gdcm::SegmentedPaletteColorLookupTable::SetLUT` ( `LookupTableType` *type*, const unsigned char \* *array*, unsigned int *length* ) [virtual]

Initialize a [SegmentedPaletteColorLookupTable](#).

Reimplemented from [gdcm::LookupTable](#).

The documentation for this class was generated from the following file:


- [gdcmSegmentedPaletteColorLookupTable.h](#)

## 27.240 `gdcm::SegmentReader` Class Reference


This class defines a segment reader. It reads attributes of group 0x0062.

```
#include <gdcmSegmentReader.h>
```

Inheritance diagram for `gdcm::SegmentReader`:

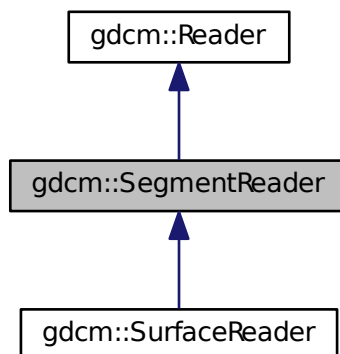

Collaboration diagram for gdcm::SegmentReader:

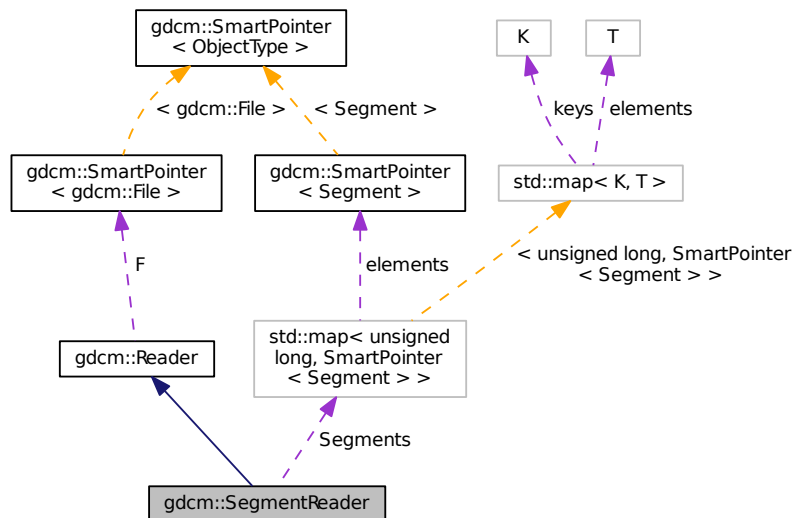

## Public Types

- typedef `std::vector< SmartPointer< Segment > >` [SegmentVector](#)

## Public Member Functions

- [SegmentReader](#) ()
  - virtual `~SegmentReader` ()
  - const [SegmentVector](#) [GetSegments](#) () const
  - [SegmentVector](#) [GetSegments](#) ()
  - virtual bool [Read](#) ()
- Read.*

## Protected Types

- typedef `std::map< unsigned long, SmartPointer< Segment > >` [SegmentMap](#)

## Protected Member Functions

- bool [ReadSegment](#) (const [Item](#) &segmentItem, const unsigned int idx)
- bool [ReadSegments](#) ()

## Protected Attributes

- [SegmentMap](#) [Segments](#)

### 27.240.1 Detailed Description

This class defines a segment reader. It reads attributes of group 0x0062.

See also

PS 3.3 C.8.20.2 and C.8.23

### 27.240.2 Member Typedef Documentation

27.240.2.1 `typedef std::map< unsigned long, SmartPointer< Segment > > gdcm::SegmentReader::SegmentMap`  
[protected]

27.240.2.2 `typedef std::vector< SmartPointer< Segment > > gdcm::SegmentReader::SegmentVector`

### 27.240.3 Constructor & Destructor Documentation

27.240.3.1 `gdcm::SegmentReader::SegmentReader ( )`

27.240.3.2 `virtual gdcm::SegmentReader::~~SegmentReader ( )` [virtual]

### 27.240.4 Member Function Documentation

27.240.4.1 `const SegmentVector gdcm::SegmentReader::GetSegments ( ) const`

27.240.4.2 `SegmentVector gdcm::SegmentReader::GetSegments ( )`

27.240.4.3 `virtual bool gdcm::SegmentReader::Read ( )` [virtual]

Read.

Reimplemented from [gdcm::Reader](#).

Reimplemented in [gdcm::SurfaceReader](#).

27.240.4.4 `bool gdcm::SegmentReader::ReadSegment ( const Item & segmentItem, const unsigned int idx )` [protected]

27.240.4.5 `bool gdcm::SegmentReader::ReadSegments ( )` [protected]

### 27.240.5 Member Data Documentation

27.240.5.1 `SegmentMap gdcm::SegmentReader::Segments` [protected]

The documentation for this class was generated from the following file:


- [gdcmSegmentReader.h](#)

## 27.241 gdcm::SegmentWriter Class Reference


This class defines a segment writer. It writes attributes of group 0x0062.

```
#include <gdcmSegmentWriter.h>
```

Inheritance diagram for gdcm::SegmentWriter:

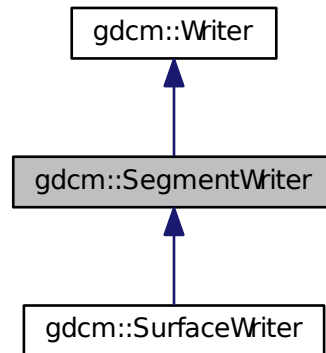

Collaboration diagram for gdcm::SegmentWriter:

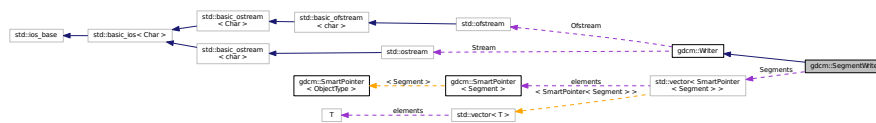

## Public Types

- typedef std::vector< [SmartPointer< Segment >](#) > [SegmentVector](#)

## Public Member Functions

- [SegmentWriter](#) ()
- virtual [~SegmentWriter](#) ()
- void [AddSegment](#) ([SmartPointer< Segment >](#) segment)
- unsigned int [GetNumberOfSegments](#) () const
- [SmartPointer< Segment >](#) [GetSegment](#) (const unsigned int idx=0) const
- const [SegmentVector](#) & [GetSegments](#) () const
- [SegmentVector](#) & [GetSegments](#) ()
- void [SetNumberOfSegments](#) (const unsigned int size)
- void [SetSegments](#) ([SegmentVector](#) &segments)
- bool [Write](#) ()

*Write.*

## Protected Member Functions

- bool [PrepareWrite](#) ()

## Protected Attributes

- [SegmentVector](#) [Segments](#)

### 27.241.1 Detailed Description

This class defines a segment writer. It writes attributes of group 0x0062.

See also

PS 3.3 C.8.20.2 and C.8.23

### 27.241.2 Member Typedef Documentation

27.241.2.1 `typedef std::vector< SmartPointer< Segment > > gdcm::SegmentWriter::SegmentVector`

### 27.241.3 Constructor & Destructor Documentation

27.241.3.1 `gdcm::SegmentWriter::SegmentWriter ( )`

27.241.3.2 `virtual gdcm::SegmentWriter::~~SegmentWriter ( )` `[virtual]`

### 27.241.4 Member Function Documentation

27.241.4.1 `void gdcm::SegmentWriter::AddSegment ( SmartPointer< Segment > segment )`

27.241.4.2 `unsigned int gdcm::SegmentWriter::GetNumberOfSegments ( )` `const`

27.241.4.3 `SmartPointer< Segment > gdcm::SegmentWriter::GetSegment ( const unsigned int idx = 0 )` `const`

27.241.4.4 `const SegmentVector& gdcm::SegmentWriter::GetSegments ( )` `const`

27.241.4.5 `SegmentVector& gdcm::SegmentWriter::GetSegments ( )`

27.241.4.6 `bool gdcm::SegmentWriter::PrepareWrite ( )` `[protected]`

27.241.4.7 `void gdcm::SegmentWriter::SetNumberOfSegments ( const unsigned int size )`

27.241.4.8 `void gdcm::SegmentWriter::SetSegments ( SegmentVector & segments )`

27.241.4.9 `bool gdcm::SegmentWriter::Write ( )` `[virtual]`

Write.

Reimplemented from [gdcm::Writer](#).

Reimplemented in [gdcm::SurfaceWriter](#).

### 27.241.5 Member Data Documentation

#### 27.241.5.1 SegmentVector gdcm::SegmentWriter::Segments [protected]

The documentation for this class was generated from the following file:


- [gdcmSegmentWriter.h](#)

## 27.242 gdcm::SequenceOfFragments Class Reference


Class to represent a Sequence Of Fragments.

```
#include <gdcmSequenceOfFragments.h>
```

Inheritance diagram for gdcm::SequenceOfFragments:

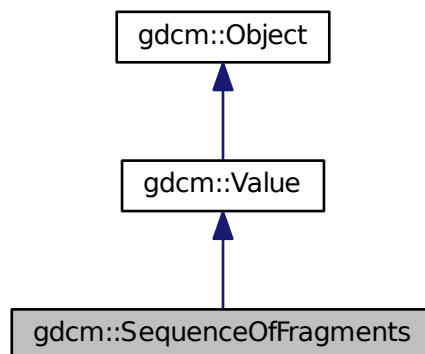

Collaboration diagram for `gdcm::SequenceOfFragments`:

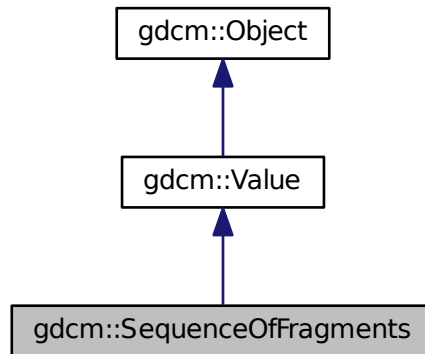

## Public Types

- typedef `FragmentVector::const_iterator` [ConstIterator](#)
- typedef `std::vector< Fragment >` [FragmentVector](#)
- typedef `FragmentVector::iterator` [Iterator](#)
- typedef `FragmentVector::size_type` [SizeType](#)

## Public Member Functions

- [SequenceOfFragments](#) ()  
*constructor (UndefinedLength by default)*
- void [AddFragment](#) ([Fragment](#) const &item)  
*Appends a [Fragment](#) to the already added ones.*
- [Iterator](#) [Begin](#) ()
- [ConstIterator](#) [Begin](#) () const
- void [Clear](#) ()  
*Clear.*
- unsigned long [ComputeByteLength](#) () const
- [VL](#) [ComputeLength](#) () const
- [Iterator](#) [End](#) ()
- [ConstIterator](#) [End](#) () const
- bool [GetBuffer](#) (char \*buffer, unsigned long length) const
- bool [GetFragBuffer](#) (unsigned int fragNb, char \*buffer, unsigned long &length) const
- const [Fragment](#) & [GetFragment](#) ([SizeType](#) num) const
- [VL](#) [GetLength](#) () const  
*Returns the SQ length, as read from disk.*
- [SizeType](#) [GetNumberOfFragments](#) () const
- const [BasicOffsetTable](#) & [GetTable](#) () const
- [BasicOffsetTable](#) & [GetTable](#) ()



- bool [operator==](#) (const [Value](#) &val) const
- void [Print](#) (std::ostream &os) const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is, bool readvalues=true)
- template<typename TSwap >  
std::istream & [ReadPreValue](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadValue](#) (std::istream &is, bool)
- void [SetLength](#) (VL length)  
*Sets the actual SQ length.*
- template<typename TSwap >  
std::ostream const & [Write](#) (std::ostream &os) const
- bool [WriteBuffer](#) (std::ostream &os) const

### Static Public Member Functions

- static [SmartPointer](#)< [SequenceOfFragments](#) > [New](#) ()

### Additional Inherited Members

#### 27.242.1 Detailed Description

Class to represent a Sequence Of Fragments.

**Todo** I do not enforce that Sequence of Fragments ends with a SQ end del

Examples:

[FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), and [GetJPEGSamplePrecision.cxx](#).

#### 27.242.2 Member Typedef Documentation

27.242.2.1 typedef [FragmentVector::const\\_iterator](#) [gdcm::SequenceOfFragments::ConstIterator](#)

27.242.2.2 typedef [std::vector<Fragment>](#) [gdcm::SequenceOfFragments::FragmentVector](#)

27.242.2.3 typedef [FragmentVector::iterator](#) [gdcm::SequenceOfFragments::Iterator](#)

27.242.2.4 typedef [FragmentVector::size\\_type](#) [gdcm::SequenceOfFragments::SizeType](#)

#### 27.242.3 Constructor & Destructor Documentation

27.242.3.1 [gdcm::SequenceOfFragments::SequenceOfFragments](#) ( ) [[inline](#)]

constructor (UndefinedLength by default)

#### 27.242.4 Member Function Documentation

27.242.4.1 void `gdcM::SequenceOfFragments::AddFragment ( Fragment const & item )`

Appends a [Fragment](#) to the already added ones.

Examples:

[FixBrokenJ2K.cxx](#).

27.242.4.2 Iterator `gdcM::SequenceOfFragments::Begin ( )` `[inline]`

27.242.4.3 ConstIterator `gdcM::SequenceOfFragments::Begin ( ) const` `[inline]`

27.242.4.4 void `gdcM::SequenceOfFragments::Clear ( )` `[virtual]`

Clear.

Implements [gdcM::Value](#).

27.242.4.5 unsigned long `gdcM::SequenceOfFragments::ComputeByteLength ( ) const`

27.242.4.6 VL `gdcM::SequenceOfFragments::ComputeLength ( ) const`

27.242.4.7 Iterator `gdcM::SequenceOfFragments::End ( )` `[inline]`

27.242.4.8 ConstIterator `gdcM::SequenceOfFragments::End ( ) const` `[inline]`

27.242.4.9 bool `gdcM::SequenceOfFragments::GetBuffer ( char * buffer, unsigned long length ) const`

27.242.4.10 bool `gdcM::SequenceOfFragments::GetFragBuffer ( unsigned int fragNb, char * buffer, unsigned long & length ) const`

27.242.4.11 const `Fragment&` `gdcM::SequenceOfFragments::GetFragment ( SizeType num ) const`

Examples:

[FixBrokenJ2K.cxx](#), and [FixJAIBugJPEGSL.cxx](#).

27.242.4.12 VL `gdcM::SequenceOfFragments::GetLength ( ) const` `[inline]`, `[virtual]`

Returns the SQ length, as read from disk.

Implements [gdcM::Value](#).

27.242.4.13 `SizeType` `gdcM::SequenceOfFragments::GetNumberOfFragments ( ) const`

Examples:

[FixJAIBugJPEGSL.cxx](#).

27.242.4.14 `const BasicOffsetTable& gdcm::SequenceOfFragments::GetTable ( ) const` `[inline]`

27.242.4.15 `BasicOffsetTable& gdcm::SequenceOfFragments::GetTable ( )` `[inline]`

27.242.4.16 `static SmartPointer<SequenceOfFragments> gdcm::SequenceOfFragments::New ( )` `[inline]`,  
`[static]`

27.242.4.17 `bool gdcm::SequenceOfFragments::operator== ( const Value & val ) const` `[inline]`,`[virtual]`

Implements [gdcm::Value](#).

27.242.4.18 `void gdcm::SequenceOfFragments::Print ( std::ostream & os ) const` `[inline]`,`[virtual]`

Reimplemented from [gdcm::Object](#).

27.242.4.19 `template<typename TSwap > std::istream& gdcm::SequenceOfFragments::Read ( std::istream & is, bool readvalues = true )` `[inline]`

27.242.4.20 `template<typename TSwap > std::istream& gdcm::SequenceOfFragments::ReadPreValue ( std::istream & is )`  
`[inline]`

References [gdcmDebugMacro](#).

27.242.4.21 `template<typename TSwap > std::istream& gdcm::SequenceOfFragments::ReadValue ( std::istream & is, bool )`  
`[inline]`

References [gdcmAssertAlwaysMacro](#), [gdcmDebugMacro](#), [gdcmWarningMacro](#), [gdcm::Tag::GetElement\(\)](#), [gdcm::Tag::GetGroup\(\)](#), [gdcm::ByteValue::GetLength\(\)](#), [gdcm::ByteValue::GetPointer\(\)](#), [gdcm::DataElement::GetTag\(\)](#), [gdcm::DataElement::GetVL\(\)](#), [gdcm::Fragment::Read\(\)](#), [gdcm::Fragment::ReadBacktrack\(\)](#), and [gdcm::Exception::what\(\)](#).

27.242.4.22 `void gdcm::SequenceOfFragments::SetLength ( VL length )` `[inline]`,`[virtual]`

Sets the actual SQ length.

Implements [gdcm::Value](#).

27.242.4.23 `template<typename TSwap > std::ostream const& gdcm::SequenceOfFragments::Write ( std::ostream & os ) const`  
`[inline]`

References [gdcm::VL::Write\(\)](#), and [gdcm::Tag::Write\(\)](#).

27.242.4.24 `bool gdcm::SequenceOfFragments::WriteBuffer ( std::ostream & os ) const`

Examples:

[GetJPEGSamplePrecision.cxx](#).

The documentation for this class was generated from the following file:


- [gdcmSequenceOfFragments.h](#)

## 27.243 gdcM::SequenceOfItems Class Reference


Class to represent a Sequence Of Items (value representation : SQ)

```
#include <gdcMSequenceOfItems.h>
```

Inheritance diagram for gdcM::SequenceOfItems:

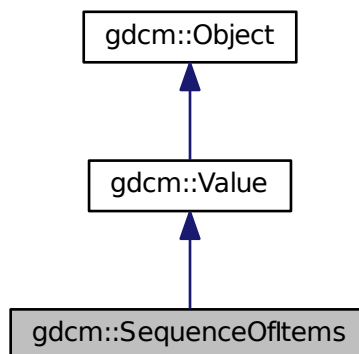

Collaboration diagram for gdcmm::SequenceOfItems:

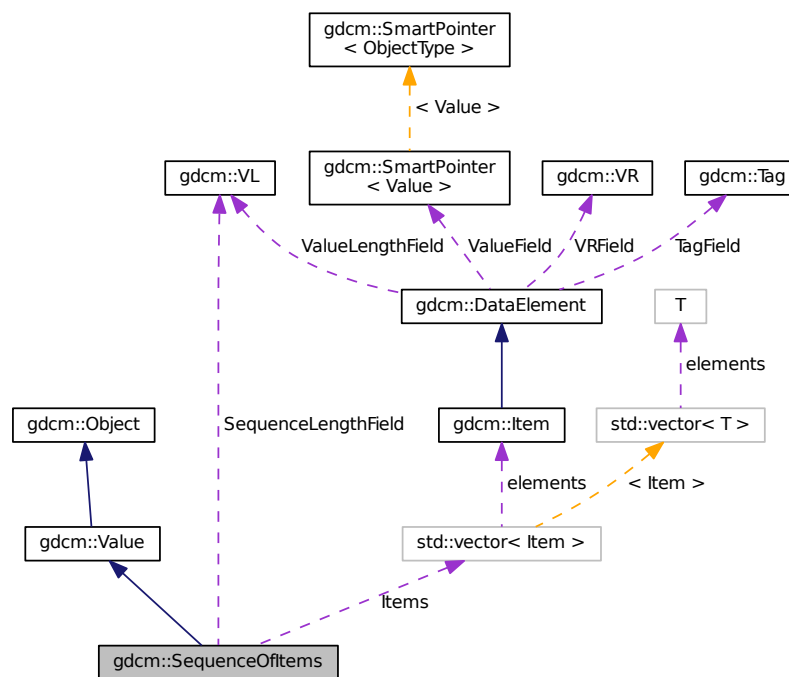

## Public Types

- `typedef ItemVector::const_iterator` [ConstIterator](#)
- `typedef std::vector< Item >` [ItemVector](#)
- `typedef ItemVector::iterator` [Iterator](#)
- `typedef ItemVector::size_type` [SizeType](#)

## Public Member Functions

- [SequenceOfItems](#) ()  
*constructor (UndefinedLength by default)*
- `void` [AddItem](#) ([Item](#) const &item)  
*Appends an Item to the already added ones.*
- [Iterator](#) [Begin](#) ()
- [ConstIterator](#) [Begin](#) () const
- `void` [Clear](#) ()  
*remove all items within the sequence*
- `template<typename TDE >`  
[VL](#) [ComputeLength](#) () const
- [Iterator](#) [End](#) ()
- [ConstIterator](#) [End](#) () const
- `bool` [FindDataElement](#) (const [Tag](#) &t) const

- const [Item](#) & [GetItem](#) ([SizeType](#) position) const
- [Item](#) & [GetItem](#) ([SizeType](#) position)
- [VL GetLength](#) () const  
*Returns the SQ length, as read from disk.*
- [SizeType GetNumberOfItems](#) () const
- bool [IsUndefinedLength](#) () const  
*return if [Value](#) Length if of undefined length*
- [SequenceOfItems](#) & [operator=](#) (const [SequenceOfItems](#) &val)
- bool [operator==](#) (const [Value](#) &val) const
- void [Print](#) (std::ostream &os) const
- template<typename TDE , typename TSwap >  
std::istream & [Read](#) (std::istream &is, bool readvalues=true)
- bool [RemoveItemByIndex](#) (const [SizeType](#) index)
- void [SetLength](#) ([VL](#) length)  
*Sets the actual SQ length.*
- void [SetLengthToUndefined](#) ()  
*Properly set the Sequence of [Item](#) to be undefined length.*
- void [SetNumberOfItems](#) ([SizeType](#) n)
- template<typename TDE , typename TSwap >  
std::ostream const & [Write](#) (std::ostream &os) const

## Static Public Member Functions

- static [SmartPointer](#)< [SequenceOfItems](#) > [New](#) ()

## Public Attributes

- [ItemVector](#) [Items](#)  
*Vector of Sequence Items.*
- [VL](#) [SequenceLengthField](#)  
*Total length of the Sequence (or 0xffffffff) if undefined.*

## Additional Inherited Members

### 27.243.1 Detailed Description

Class to represent a Sequence Of Items (value representation : SQ)

- a [Value](#) Representation for Data Elements that contains a sequence of Data Sets.
- Sequence of [Item](#) allows for Nested Data Sets

See PS 3.5, 7.4.6 Data [Element Type](#) Within a Sequence

## Note

SEQUENCE OF ITEMS (VALUE REPRESENTATION SQ) A [Value](#) Representation for Data Elements that contain a sequence of Data Sets. Sequence of Items allows for Nested Data Sets.

## Examples:

[DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSequenceUltrasound.cxx](#), and [ReadExplicitLengthSQIVR.cxx](#).

## 27.243.2 Member Typedef Documentation

27.243.2.1 `typedef ItemVector::const_iterator gdcm::SequenceOfItems::ConstIterator`

27.243.2.2 `typedef std::vector< Item > gdcm::SequenceOfItems::ItemVector`

27.243.2.3 `typedef ItemVector::iterator gdcm::SequenceOfItems::Iterator`

27.243.2.4 `typedef ItemVector::size_type gdcm::SequenceOfItems::SizeType`

## 27.243.3 Constructor &amp; Destructor Documentation

27.243.3.1 `gdcm::SequenceOfItems::SequenceOfItems ( ) [inline]`

constructor (UndefinedLength by default)

## 27.243.4 Member Function Documentation

27.243.4.1 `void gdcm::SequenceOfItems::AddItem ( Item const & item )`

Appends an [Item](#) to the already added ones.

## Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenLongSeqs.cxx](#), and [GenSeqs.cxx](#).

27.243.4.2 `Iterator gdcm::SequenceOfItems::Begin ( ) [inline]`

27.243.4.3 `ConstIterator gdcm::SequenceOfItems::Begin ( ) const [inline]`

27.243.4.4 `void gdcm::SequenceOfItems::Clear ( ) [virtual]`

remove all items within the sequence

Implements [gdcm::Value](#).

27.243.4.5 `template<typename TDE > VL gdcm::SequenceOfItems::ComputeLength ( ) const`

27.243.4.6 `Iterator gdcm::SequenceOfItems::End ( ) [inline]`

27.243.4.7 **ConstIterator** `gdcm::SequenceOfItems::End ( ) const` `[inline]`

27.243.4.8 **bool** `gdcm::SequenceOfItems::FindDataElement ( const Tag & t ) const`

27.243.4.9 **const Item&** `gdcm::SequenceOfItems::GetItem ( SizeType position ) const`

Examples:

[ChangeSequenceUltrasound.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [ExtractEncryptedContent.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [LargeVRDSExplicit.cxx](#), and [ReadAndDumpDICOMDIR.cxx](#).

27.243.4.10 **Item&** `gdcm::SequenceOfItems::GetItem ( SizeType position )`

27.243.4.11 **VL** `gdcm::SequenceOfItems::GetLength ( ) const` `[inline],[virtual]`

Returns the SQ length, as read from disk.

Implements [gdcm::Value](#).

27.243.4.12 **SizeType** `gdcm::SequenceOfItems::GetNumberOfItems ( ) const` `[inline]`

Examples:

[ChangeSequenceUltrasound.cxx](#), [DumpExamCard.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECHO.cxx](#), [ExtractEncryptedContent.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), and [LargeVRDSExplicit.cxx](#).

27.243.4.13 **bool** `gdcm::SequenceOfItems::IsUndefinedLength ( ) const` `[inline]`

return if [Value](#) Length if of undefined length

27.243.4.14 **static SmartPointer<SequenceOfItems>** `gdcm::SequenceOfItems::New ( )` `[inline],[static]`

Examples:

[NewSequence.cs](#).

27.243.4.15 **SequenceOfItems&** `gdcm::SequenceOfItems::operator= ( const SequenceOfItems & val )` `[inline]`

References Items, and SequenceLengthField.

27.243.4.16 **bool** `gdcm::SequenceOfItems::operator==( const Value & val ) const` `[inline],[virtual]`

Implements [gdcm::Value](#).

References Items, and SequenceLengthField.



27.243.4.17 `void gdcm::SequenceOfItems::Print ( std::ostream & os ) const` `[inline]`, `[virtual]`

Reimplemented from [gdcm::Object](#).

27.243.4.18 `template<typename TDE , typename TSwap > std::istream& gdcm::SequenceOfItems::Read ( std::istream & is, bool readvalues = true )` `[inline]`

Examples:

[ReadExplicitLengthSQIVR.cxx](#).

References [gdcm::Item::Clear\(\)](#), [gdcmDebugMacro](#), [gdcmWarningMacro](#), [gdcm::Exception::GetDescription\(\)](#), [gdcm::Item::GetNestedDataSet\(\)](#), [gdcm::DataElement::GetTag\(\)](#), [gdcm::DataElement::GetVL\(\)](#), [gdcm::Item::Read\(\)](#), and [gdcm::DataSet::Size\(\)](#).

27.243.4.19 `bool gdcm::SequenceOfItems::RemoveItemByIndex ( const SizeType index )`

Remove an [Item](#) as specified by its index, if index > size, false is returned Index starts at 1 not 0

27.243.4.20 `void gdcm::SequenceOfItems::SetLength ( VL length )` `[inline]`, `[virtual]`

Sets the actual SQ length.

Implements [gdcm::Value](#).

Examples:

[ReadExplicitLengthSQIVR.cxx](#).

27.243.4.21 `void gdcm::SequenceOfItems::SetLengthToUndefined ( )`

Properly set the Sequence of [Item](#) to be undefined length.

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllIVR.cxx](#), [GenLongSeqs.cxx](#), and [GenSeqs.cxx](#).

27.243.4.22 `void gdcm::SequenceOfItems::SetNumberOfItems ( SizeType n )` `[inline]`

27.243.4.23 `template<typename TDE , typename TSwap > std::ostream const& gdcm::SequenceOfItems::Write ( std::ostream & os ) const` `[inline]`

References [gdcm::VL::Write\(\)](#), and [gdcm::Tag::Write\(\)](#).

## 27.243.5 Member Data Documentation

### 27.243.5.1 ItemVector [gdcm::SequenceOfItems::Items](#)

Vector of Sequence Items.

Referenced by [operator=\(\)](#), and [operator==\(\)](#).

### 27.243.5.2 VL gdcm::SequenceOfItems::SequenceLengthField

Total length of the Sequence (or 0xffffffff if undefined).

Referenced by operator=(), and operator==().

The documentation for this class was generated from the following file:


- [gdcmSequenceOfItems.h](#)

## 27.244 gdcm::SerieHelper Class Reference

[SerieHelper](#) DO NOT USE this class, it is only a temporary solution for ITK migration from GDCM 1.x to GDCM 2.x It will disappear soon, you've been warned.

```
#include <gdcmSerieHelper.h>
```

Collaboration diagram for gdcm::SerieHelper:

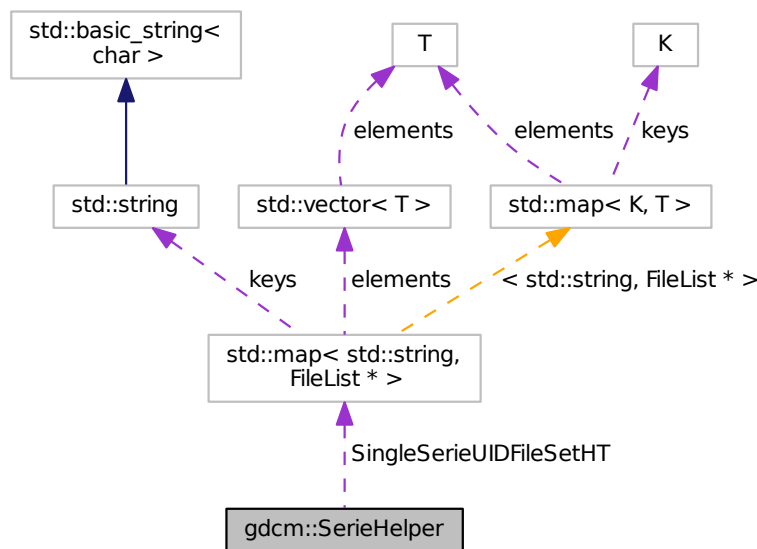

### Classes

- struct [Rule](#)

### Public Member Functions

- [SerieHelper](#) ()
- [~SerieHelper](#) ()
- void [AddRestriction](#) (const std::string &tag)

- void [AddRestriction](#) (uint16\_t group, uint16\_t elem, std::string const &value, int op)
- void [Clear](#) ()
- void [CreateDefaultUniqueSeriesIdentifier](#) ()
- std::string [CreateUniqueSeriesIdentifier](#) ([File](#) \*inFile)
- [FileList](#) \* [GetFirstSingleSerieUIDFileSet](#) ()
- [FileList](#) \* [GetNextSingleSerieUIDFileSet](#) ()
- void [OrderFileList](#) ([FileList](#) \*fileSet)
- void [SetDirectory](#) (std::string const &dir, bool recursive=false)
- void [SetLoadMode](#) (int)
- void [SetUseSeriesDetails](#) (bool useSeriesDetails)

## Protected Types

- typedef std::vector< [Rule](#) > [SerieRestrictions](#)
- typedef std::map< std::string, [FileList](#) \* > [SingleSerieUIDFileSetmap](#)

## Protected Member Functions

- bool [AddFile](#) ([FileWithName](#) &header)
- void [AddFileName](#) (std::string const &filename)
- void [AddRestriction](#) (const [Tag](#) &tag)
- bool [FileNameOrdering](#) ([FileList](#) \*fileList)
- bool [ImagePositionPatientOrdering](#) ([FileList](#) \*fileSet)
- bool [UserOrdering](#) ([FileList](#) \*fileSet)

## Protected Attributes

- [SingleSerieUIDFileSetmap::iterator](#) [ItFileSetHt](#)
- [SingleSerieUIDFileSetmap](#) [SingleSerieUIDFileSetHT](#)

### 27.244.1 Detailed Description

[SerieHelper](#) DO NOT USE this class, it is only a temporary solution for ITK migration from GDCM 1.x to GDCM 2.x It will disappear soon, you've been warned.

Instead see [ImageHelper](#) or [IPPSorter](#)

### 27.244.2 Member Typedef Documentation

27.244.2.1 typedef std::vector<[Rule](#)> [gdcm::SerieHelper::SerieRestrictions](#) [protected]

27.244.2.2 typedef std::map<std::string, [FileList](#) \*> [gdcm::SerieHelper::SingleSerieUIDFileSetmap](#) [protected]

### 27.244.3 Constructor & Destructor Documentation

27.244.3.1 [gdcm::SerieHelper::SerieHelper](#) ( )

27.244.3.2 `gdcm::SerieHelper::~~SerieHelper ( )`

## 27.244.4 Member Function Documentation

27.244.4.1 `bool gdcm::SerieHelper::AddFile ( FileWithName & header )` [protected]

27.244.4.2 `void gdcm::SerieHelper::AddFileName ( std::string const & filename )` [protected]

27.244.4.3 `void gdcm::SerieHelper::AddRestriction ( const std::string & tag )`

27.244.4.4 `void gdcm::SerieHelper::AddRestriction ( uint16_t group, uint16_t elem, std::string const & value, int op )`

27.244.4.5 `void gdcm::SerieHelper::AddRestriction ( const Tag & tag )` [protected]

27.244.4.6 `void gdcm::SerieHelper::Clear ( )`

27.244.4.7 `void gdcm::SerieHelper::CreateDefaultUniqueSeriesIdentifier ( )`

27.244.4.8 `std::string gdcm::SerieHelper::CreateUniqueSeriesIdentifier ( File * inFile )`

27.244.4.9 `bool gdcm::SerieHelper::FileNameOrdering ( FileList * fileList )` [protected]

27.244.4.10 `FileList* gdcm::SerieHelper::GetFirstSingleSerieUIDFileSet ( )`

27.244.4.11 `FileList* gdcm::SerieHelper::GetNextSingleSerieUIDFileSet ( )`

27.244.4.12 `bool gdcm::SerieHelper::ImagePositionPatientOrdering ( FileList * fileSet )` [protected]

27.244.4.13 `void gdcm::SerieHelper::OrderFileList ( FileList * fileSet )`

27.244.4.14 `void gdcm::SerieHelper::SetDirectory ( std::string const & dir, bool recursive = false )`

27.244.4.15 `void gdcm::SerieHelper::SetLoadMode ( int )` [inline]

27.244.4.16 `void gdcm::SerieHelper::SetUseSeriesDetails ( bool useSeriesDetails )`

27.244.4.17 `bool gdcm::SerieHelper::UserOrdering ( FileList * fileSet )` [protected]

## 27.244.5 Member Data Documentation

27.244.5.1 `SingleSerieUIDFileSetmap::iterator gdcm::SerieHelper::ItFileSetHt` [protected]

27.244.5.2 `SingleSerieUIDFileSetmap gdcm::SerieHelper::SingleSerieUIDFileSetHT` [protected]

The documentation for this class was generated from the following file:

- [gdcmSerieHelper.h](#)

## 27.245 gdcm::Series Class Reference

[Series.](#)

```
#include <gdcmSeries.h>
```

## Public Member Functions

- [Series](#) ()

### 27.245.1 Detailed Description

[Series](#).

### 27.245.2 Constructor & Destructor Documentation

27.245.2.1 `gdcm::Series::Series ( )` `[inline]`

The documentation for this class was generated from the following file:

- [gdcmSeries.h](#)

## 27.246 gdcm::network::ServiceClassApplicationInformation Class Reference

```
#include <gdcmServiceClassApplicationInformation.h>
```

## Public Member Functions

- [ServiceClassApplicationInformation](#) ()
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetTuple](#) (uint8\_t levelofsupport, uint8\_t levelofdigitalsig, uint8\_t elementcoercion)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.246.1 Detailed Description

PS 3.4 [Table B.3-1](#) SERVICE-CLASS-APPLICATION-INFORMATION (A-ASSOCIATE-RQ)

### 27.246.2 Constructor & Destructor Documentation

27.246.2.1 `gdcm::network::ServiceClassApplicationInformation::ServiceClassApplicationInformation ( )`

### 27.246.3 Member Function Documentation

27.246.3.1 `void gdcm::network::ServiceClassApplicationInformation::Print ( std::ostream & os ) const`

27.246.3.2 `std::istream& gdcm::network::ServiceClassApplicationInformation::Read ( std::istream & is )`

27.246.3.3 void gdcmm::network::ServiceClassApplicationInformation::SetTuple ( uint8\_t *levelofsupport*, uint8\_t *levelofdigitalsig*,  
uint8\_t *elementcoercion* )

27.246.3.4 size\_t gdcmm::network::ServiceClassApplicationInformation::Size ( ) const

27.246.3.5 const std::ostream& gdcmm::network::ServiceClassApplicationInformation::Write ( std::ostream & *os* ) const

The documentation for this class was generated from the following file:


- [gdcmmServiceClassApplicationInformation.h](#)

## 27.247 gdcmm::ServiceClassUser Class Reference


[ServiceClassUser](#).

```
#include <gdcmmServiceClassUser.h>
```

Inheritance diagram for gdcmm::ServiceClassUser:

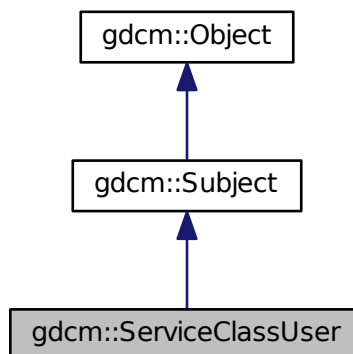

Collaboration diagram for gdcm::ServiceClassUser:

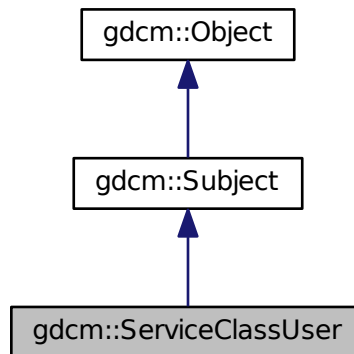

## Public Member Functions

- [ServiceClassUser](#) ()
- [~ServiceClassUser](#) ()
- const char \* [GetAETitle](#) () const
- const char \* [GetCalledAETitle](#) () const
- double [GetTimeout](#) () const
- bool [InitializeConnection](#) ()
- bool [IsPresentationContextAccepted](#) (const [PresentationContext](#) &pc) const  
*Return if the passed in presentation was accepted during association negotiation.*
- bool [SendEcho](#) ()  
*C-ECHO.*
- bool [SendFind](#) (const [BaseRootQuery](#) \*query, std::vector< [DataSet](#) > &retDatasets)  
*C-FIND a query, return result are in retDatasets.*
- bool [SendMove](#) (const [BaseRootQuery](#) \*query, const char \*outputdir)  
*Execute a C-MOVE, based on query, return files are written in outputdir.*
- bool [SendMove](#) (const [BaseRootQuery](#) \*query, std::vector< [DataSet](#) > &retDatasets)  
*Execute a C-MOVE, based on query, returned dataset are Implicit.*
- bool [SendMove](#) (const [BaseRootQuery](#) \*query, std::vector< [File](#) > &retFile)  
*Execute a C-MOVE, based on query, returned Files are stored in vector.*
- bool [SendStore](#) (const char \*filename)  
*Execute a C-STORE on file on disk, named filename.*
- bool [SendStore](#) ([File](#) const &file)
- bool [SendStore](#) ([DataSet](#) const &ds)  
*Execute a C-STORE on a DataSet, the transfer syntax used will be Implicit.*
- void [SetAETitle](#) (const char \*aetitle)  
*set calling ae title*
- void [SetCalledAETitle](#) (const char \*aetitle)

- set called ae title*
- void [SetHostname](#) (const char \*hostname)
  - Set the name of the called hostname (hostname or IP address)*
- void [SetPort](#) (uint16\_t port)
  - Set port of remote host (called application)*
- void [SetPortSCP](#) (uint16\_t portscp)
  - Set the port for any incoming C-STORE-SCP operation (typically in a return of C-MOVE)*
- void [SetPresentationContexts](#) (std::vector< [PresentationContext](#) > const &pcs)
  - Set the Presentation Context used for the Association.*
- void [SetTimeout](#) (double t)
  - set/get Timeout*
- bool [StartAssociation](#) ()
  - Start the association. Need to call SetPresentationContexts before.*
- bool [StopAssociation](#) ()
  - Stop the running association.*

## Static Public Member Functions

- static [SmartPointer](#)< [ServiceClassUser](#) > [New](#) ()
  - for wrapped language: instantiate a reference counted object*

## Additional Inherited Members

### 27.247.1 Detailed Description

[ServiceClassUser](#).

Examples:

[CStoreQtProgress.cxx](#).

### 27.247.2 Constructor & Destructor Documentation

#### 27.247.2.1 `gdcm::ServiceClassUser::ServiceClassUser ( )`

Construct a SCU with default:

- hostname = localhost
- port = 104

#### 27.247.2.2 `gdcm::ServiceClassUser::~~ServiceClassUser ( )`

### 27.247.3 Member Function Documentation

#### 27.247.3.1 `const char* gdcm::ServiceClassUser::GetAETitle ( ) const`



27.247.3.2 `const char* gdcm::ServiceClassUser::GetCalledAETitle ( ) const`

27.247.3.3 `double gdcm::ServiceClassUser::GetTimeout ( ) const`

27.247.3.4 `bool gdcm::ServiceClassUser::InitializeConnection ( )`

Will try to connect This will setup the actual timeout used during the whole connection time. Need to call SetTimeout first

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.5 `bool gdcm::ServiceClassUser::IsPresentationContextAccepted ( const PresentationContext & pc ) const`

Return if the passed in presentation was accepted during association negotiation.

27.247.3.6 `static SmartPointer<ServiceClassUser> gdcm::ServiceClassUser::New ( ) [inline],[static]`

for wrapped language: instantiate a reference counted object

27.247.3.7 `bool gdcm::ServiceClassUser::SendEcho ( )`

C-ECHO.

27.247.3.8 `bool gdcm::ServiceClassUser::SendFind ( const BaseRootQuery * query, std::vector< DataSet > & retDatasets )`

C-FIND a query, return result are in retDatasets.

27.247.3.9 `bool gdcm::ServiceClassUser::SendMove ( const BaseRootQuery * query, const char * outputdir )`

Execute a C-MOVE, based on query, return files are written in outputdir.

27.247.3.10 `bool gdcm::ServiceClassUser::SendMove ( const BaseRootQuery * query, std::vector< DataSet > & retDatasets )`

Execute a C-MOVE, based on query, returned dataset are Implicit.

27.247.3.11 `bool gdcm::ServiceClassUser::SendMove ( const BaseRootQuery * query, std::vector< File > & retFile )`

Execute a C-MOVE, based on query, returned Files are stored in vector.

27.247.3.12 `bool gdcm::ServiceClassUser::SendStore ( const char * filename )`

Execute a C-STORE on file on disk, named filename.

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.13 `bool gdcmm::ServiceClassUser::SendStore ( File const & file )`

Execute a C-STORE on a [File](#), the transfer syntax used for the query is based on the file.

27.247.3.14 `bool gdcmm::ServiceClassUser::SendStore ( DataSet const & ds )`

Execute a C-STORE on a [DataSet](#), the transfer syntax used will be Implicit.

27.247.3.15 `void gdcmm::ServiceClassUser::SetAETitle ( const char * aetitle )`

set calling ae title

27.247.3.16 `void gdcmm::ServiceClassUser::SetCalledAETitle ( const char * aetitle )`

set called ae title

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.17 `void gdcmm::ServiceClassUser::SetHostname ( const char * hostname )`

Set the name of the called hostname (hostname or IP address)

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.18 `void gdcmm::ServiceClassUser::SetPort ( uint16_t port )`

Set port of remote host (called application)

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.19 `void gdcmm::ServiceClassUser::SetPortSCP ( uint16_t portscp )`

Set the port for any incoming C-STORE-SCP operation (typically in a return of C-MOVE)

27.247.3.20 `void gdcmm::ServiceClassUser::SetPresentationContexts ( std::vector< PresentationContext > const & pcs )`

Set the Presentation Context used for the Association.

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.21 void gdcm::ServiceClassUser::SetTimeout ( double t )

set/get Timeout

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.22 bool gdcm::ServiceClassUser::StartAssociation ( )

Start the association. Need to call SetPresentationContexts before.

Examples:

[CStoreQtProgress.cxx](#).

27.247.3.23 bool gdcm::ServiceClassUser::StopAssociation ( )

Stop the running association.

Examples:

[CStoreQtProgress.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmServiceClassUser.h](#)

## 27.248 gdcm::SHA1 Class Reference

Class for [SHA1](#).

```
#include <gdcmSHA1.h>
```

### Public Member Functions

- [SHA1](#) ()
- [~SHA1](#) ()

### Static Public Member Functions

- static bool [Compute](#) (const char \*buffer, unsigned long buf\_len, char digest\_str[20 \*2+1])
- static bool [ComputeFile](#) (const char \*filename, char digest\_str[20 \*2+1])

### 27.248.1 Detailed Description

Class for [SHA1](#).

#### Warning

this class is able to pick from one implementation:

1. the one from OpenSSL (when GDCM\_USE\_SYSTEM\_OPENSSL is turned ON)

In all other cases it will return an error

### 27.248.2 Constructor & Destructor Documentation

27.248.2.1 `gdcm::SHA1::SHA1 ( )`

27.248.2.2 `gdcm::SHA1::~~SHA1 ( )`

### 27.248.3 Member Function Documentation

27.248.3.1 `static bool gdcm::SHA1::Compute ( const char * buffer, unsigned long buf_len, char digest_str[20 * 2 + 1] )`  
[static]

27.248.3.2 `static bool gdcm::SHA1::ComputeFile ( const char * filename, char digest_str[20 * 2 + 1] )` [static]

The documentation for this class was generated from the following file:


- [gdcmSHA1.h](#)

## 27.249 `gdcm::SimpleMemberCommand< T >` Class Template Reference


[Command](#) subclass that calls a pointer to a member function.

```
#include <gdcmCommand.h>
```

Inheritance diagram for gdcmm::SimpleMemberCommand< T >:

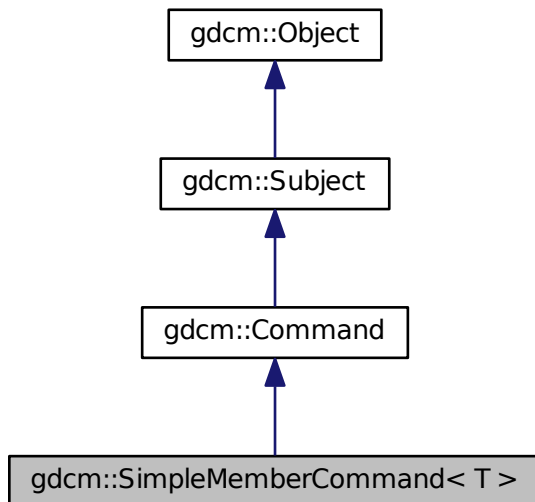

Collaboration diagram for gdcmm::SimpleMemberCommand< T >:

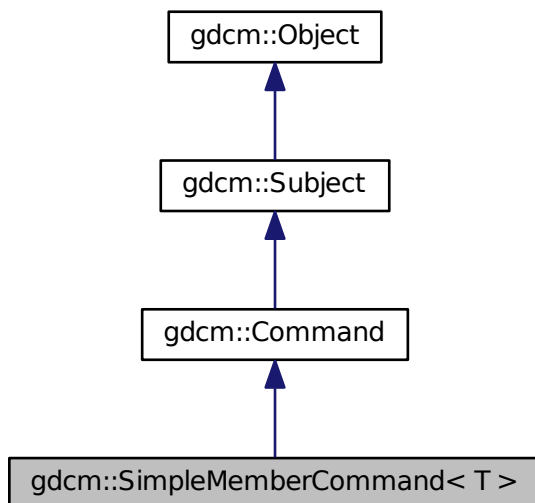

## Public Types

- typedef [SimpleMemberCommand](#) Self
- typedef void(T::\* [TMemberFunctionPointer](#)) ()

## Public Member Functions

- virtual void [Execute](#) ([Subject](#) \*, const [Event](#) &)
- virtual void [Execute](#) (const [Subject](#) \*, const [Event](#) &)
- void [SetCallbackFunction](#) (T \*object, [TMemberFunctionPointer](#) memberFunction)

## Static Public Member Functions

- static [SmartPointer](#)< [SimpleMemberCommand](#) > [New](#) ()

## Protected Member Functions

- [SimpleMemberCommand](#) ()
- virtual [~SimpleMemberCommand](#) ()

## Protected Attributes

- [TMemberFunctionPointer](#) m\_MemberFunction
- T \* m\_This

### 27.249.1 Detailed Description

template<typename T>class gdcM::SimpleMemberCommand< T >

[Command](#) subclass that calls a pointer to a member function.

[SimpleMemberCommand](#) calls a pointer to a member function with no arguments.

### 27.249.2 Member Typedef Documentation

27.249.2.1 template<typename T > typedef [SimpleMemberCommand](#) gdcM::SimpleMemberCommand< T >::Self

Standard class typedefs.

27.249.2.2 template<typename T > typedef void(T::\* [gdcM::SimpleMemberCommand](#)< T >::TMemberFunctionPointer) ()

A method callback.

### 27.249.3 Constructor & Destructor Documentation

27.249.3.1 `template<typename T> gdcmm::SimpleMemberCommand< T >::SimpleMemberCommand ( )`  
`[inline], [protected]`

Referenced by `gdcmm::SimpleMemberCommand< T >::New()`.

27.249.3.2 `template<typename T> virtual gdcmm::SimpleMemberCommand< T >::~~SimpleMemberCommand ( )`  
`[inline], [protected], [virtual]`

### 27.249.4 Member Function Documentation

27.249.4.1 `template<typename T> virtual void gdcmm::SimpleMemberCommand< T >::Execute ( Subject *, const Event & )` `[inline], [virtual]`

Invoke the callback function.

Implements [gdcmm::Command](#).

References `gdcmm::SimpleMemberCommand< T >::m_MemberFunction`.

27.249.4.2 `template<typename T> virtual void gdcmm::SimpleMemberCommand< T >::Execute ( const Subject * caller, const Event & event )` `[inline], [virtual]`

Abstract method that defines the action to be taken by the command. This variant is expected to be used when requests comes from a const [Object](#)

Implements [gdcmm::Command](#).

References `gdcmm::SimpleMemberCommand< T >::m_MemberFunction`.

27.249.4.3 `template<typename T> static SmartPointer<SimpleMemberCommand> gdcmm::SimpleMemberCommand< T >::New ( )` `[inline], [static]`

Run-time type information (and related methods). Method for creation through the object factory.

References `gdcmm::SimpleMemberCommand< T >::SimpleMemberCommand()`.

27.249.4.4 `template<typename T> void gdcmm::SimpleMemberCommand< T >::SetCallbackFunction ( T * object, TMemberFunctionPointer memberFunction )` `[inline]`

Specify the callback function.

References `gdcmm::SimpleMemberCommand< T >::m_MemberFunction`, and `gdcmm::SimpleMemberCommand< T >::m_This`.

### 27.249.5 Member Data Documentation

27.249.5.1 `template<typename T> TMemberFunctionPointer gdcmm::SimpleMemberCommand< T >::m_MemberFunction` `[protected]`

Referenced by `gdcmm::SimpleMemberCommand< T >::Execute()`, and `gdcmm::SimpleMemberCommand< T >::SetCallbackFunction()`.

27.249.5.2 `template<typename T> T* gdcm::SimpleMemberCommand< T >::m_This` [protected]

Referenced by `gdcm::SimpleMemberCommand< T >::SetCallbackFunction()`.

The documentation for this class was generated from the following file:

- [gdcmCommand.h](#)

## 27.250 gdcm::SimpleSubjectWatcher Class Reference

[SimpleSubjectWatcher](#) This is a typical [Subject](#) Watcher class. It will observe all events.

```
#include <gdcmSimpleSubjectWatcher.h>
```

### Public Member Functions

- [SimpleSubjectWatcher](#) ([Subject](#) \*s, const char \*comment="")
- virtual [~SimpleSubjectWatcher](#) ()

### Protected Member Functions

- virtual void [EndFilter](#) ()
- virtual void [ShowAbort](#) ()
- virtual void [ShowAnonymization](#) ([Subject](#) \*caller, const [Event](#) &evt)
- virtual void [ShowData](#) ([Subject](#) \*caller, const [Event](#) &evt)
- virtual void [ShowDataSet](#) ([Subject](#) \*caller, const [Event](#) &evt)
- virtual void [ShowFileName](#) ([Subject](#) \*caller, const [Event](#) &evt)
- virtual void [ShowIteration](#) ()
- virtual void [ShowProgress](#) ([Subject](#) \*caller, const [Event](#) &evt)
- virtual void [StartFilter](#) ()
- void [TestAbortOff](#) ()
- void [TestAbortOn](#) ()

### 27.250.1 Detailed Description

[SimpleSubjectWatcher](#) This is a typical [Subject](#) Watcher class. It will observe all events.

Examples:

[SimpleScanner.cxx](#).

### 27.250.2 Constructor & Destructor Documentation

27.250.2.1 `gdcm::SimpleSubjectWatcher::SimpleSubjectWatcher ( Subject * s, const char * comment = " " )`

27.250.2.2 `virtual gdcm::SimpleSubjectWatcher::~~SimpleSubjectWatcher ( )` [virtual]

### 27.250.3 Member Function Documentation



- 27.250.3.1 virtual void gdcm::SimpleSubjectWatcher::EndFilter ( ) [protected],[virtual]
- 27.250.3.2 virtual void gdcm::SimpleSubjectWatcher::ShowAbort ( ) [protected],[virtual]
- 27.250.3.3 virtual void gdcm::SimpleSubjectWatcher::ShowAnonymization ( Subject \* *caller*, const Event & *evt* ) [protected],[virtual]
- 27.250.3.4 virtual void gdcm::SimpleSubjectWatcher::ShowData ( Subject \* *caller*, const Event & *evt* ) [protected],[virtual]
- 27.250.3.5 virtual void gdcm::SimpleSubjectWatcher::ShowDataSet ( Subject \* *caller*, const Event & *evt* ) [protected],[virtual]
- 27.250.3.6 virtual void gdcm::SimpleSubjectWatcher::ShowFileName ( Subject \* *caller*, const Event & *evt* ) [protected],[virtual]

Examples:

[SimpleScanner.cxx](#).

- 27.250.3.7 virtual void gdcm::SimpleSubjectWatcher::ShowIteration ( ) [protected],[virtual]
- 27.250.3.8 virtual void gdcm::SimpleSubjectWatcher::ShowProgress ( Subject \* *caller*, const Event & *evt* ) [protected],[virtual]
- 27.250.3.9 virtual void gdcm::SimpleSubjectWatcher::StartFilter ( ) [protected],[virtual]
- 27.250.3.10 void gdcm::SimpleSubjectWatcher::TestAbortOff ( ) [protected]
- 27.250.3.11 void gdcm::SimpleSubjectWatcher::TestAbortOn ( ) [protected]

The documentation for this class was generated from the following file:


- [gdcmSimpleSubjectWatcher.h](#)

## 27.251 gdcm::SmartPointer< ObjectType > Class Template Reference

Class for Smart Pointer.

```
#include <gdcmObject.h>
```

Inheritance diagram for `gdcmm::SmartPointer< ObjectType >`:

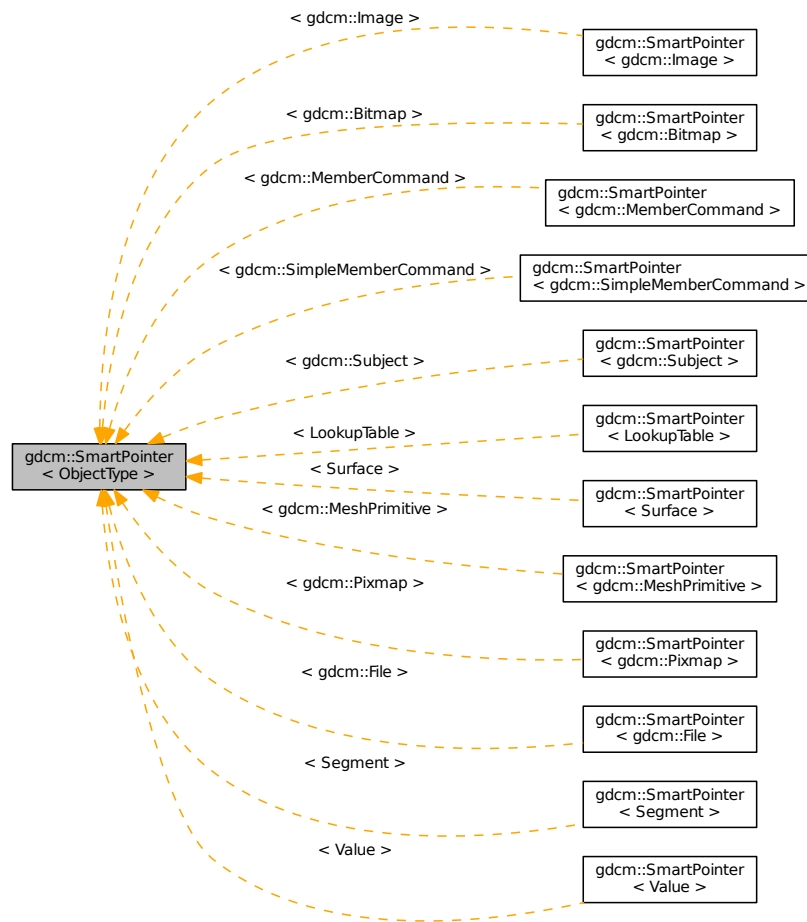

## Public Member Functions

- `SmartPointer ()`
- `SmartPointer (const SmartPointer< ObjectType > &p)`
- `SmartPointer (ObjectType *p)`
- `SmartPointer (ObjectType const &p)`
- `~SmartPointer ()`
- `ObjectType * GetPointer () const`  
*Explicit function to retrieve the pointer.*
- `operator ObjectType * () const`  
*Return pointer to object.*
- `ObjectType & operator* () const`
- `ObjectType * operator-> () const`  
*Overload operator ->*
- `SmartPointer & operator= (SmartPointer const &r)`

*Overload operator assignment.*

- [SmartPointer](#) & [operator=](#) (ObjectType \*r)

*Overload operator assignment.*

- [SmartPointer](#) & [operator=](#) (ObjectType const &r)

### 27.251.1 Detailed Description

```
template<class ObjectType>class gdcm::SmartPointer< ObjectType >
```

Class for Smart Pointer.

Will only work for subclass of [gdcm::Object](#) See `tr1/shared_ptr` for a more general approach (not invasive) `#include <tr1/memory> { shared_ptr<Bla> b(new Bla); }`

#### Note

Class partly based on post by Bill Hubauer: <http://groups.google.com/group/comp.lang.c++.msg/173ddc38a827a930>

#### See also

<http://www.davethehat.com/articles/smartp.htm>

and `itk::SmartPointer`

#### Examples:

[ChangeSequenceUltrasound.cxx](#), [CStoreQtProgress.cxx](#), [DumpGEMSMovieGroup.cxx](#), [DumpPhilipsECGHO.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSubSequenceData.cxx](#), [LargeVRDSExplicit.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), and [SimpleScanner.cxx](#).

### 27.251.2 Constructor & Destructor Documentation

27.251.2.1 `template<class ObjectType> gdcm::SmartPointer< ObjectType >::SmartPointer ( ) [inline]`

27.251.2.2 `template<class ObjectType> gdcm::SmartPointer< ObjectType >::SmartPointer ( const SmartPointer< ObjectType > & p ) [inline]`

27.251.2.3 `template<class ObjectType> gdcm::SmartPointer< ObjectType >::SmartPointer ( ObjectType * p ) [inline]`

27.251.2.4 `template<class ObjectType> gdcm::SmartPointer< ObjectType >::SmartPointer ( ObjectType const & p ) [inline]`

27.251.2.5 `template<class ObjectType> gdcm::SmartPointer< ObjectType >::~~SmartPointer ( ) [inline]`

### 27.251.3 Member Function Documentation

27.251.3.1 `template<class ObjectType> ObjectType* gdcm::SmartPointer< ObjectType >::GetPointer ( ) const [inline]`

Explicit function to retrieve the pointer.

27.251.3.2 `template<class ObjectType> gdcmm::SmartPointer< ObjectType >::operator ObjectType * ( ) const`  
`[inline]`

Return pointer to object.

27.251.3.3 `template<class ObjectType> ObjectType& gdcmm::SmartPointer< ObjectType >::operator* ( ) const`  
`[inline]`

27.251.3.4 `template<class ObjectType> ObjectType* gdcmm::SmartPointer< ObjectType >::operator-> ( ) const`  
`[inline]`

Overload operator ->

27.251.3.5 `template<class ObjectType> SmartPointer& gdcmm::SmartPointer< ObjectType >::operator= ( SmartPointer< ObjectType > const & r )` `[inline]`

Overload operator assignment.

Referenced by `gdcmm::SmartPointer< Value >::operator=()`.

27.251.3.6 `template<class ObjectType> SmartPointer& gdcmm::SmartPointer< ObjectType >::operator= ( ObjectType * r )`  
`[inline]`

Overload operator assignment.

27.251.3.7 `template<class ObjectType> SmartPointer& gdcmm::SmartPointer< ObjectType >::operator= ( ObjectType const & r )` `[inline]`

The documentation for this class was generated from the following files:

- [gdcmmObject.h](#)
- [gdcmmSmartPointer.h](#)

## 27.252 gdcmm::network::SOPClassExtendedNegociationSub Class Reference

[SOPClassExtendedNegociationSub](#) PS 3.7 [Table D.3-11](#) SOP CLASS EXTENDED NEGOTIATION SUB-ITEM FIELDS (A-ASSOCIATE-RQ and A-ASSOCIATE-AC)

```
#include <gdcmmSOPClassExtendedNegociationSub.h>
```

### Public Member Functions

- [SOPClassExtendedNegociationSub](#) ()
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- void [SetTuple](#) (const char \*uid, uint8\_t levelofsupport=3, uint8\_t levelofdignalsig=0, uint8\_t elementcoercion=2)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.252.1 Detailed Description

[SOPClassExtendedNegociationSub](#) PS 3.7 [Table D.3-11](#) SOP CLASS EXTENDED NEGOTIATION SUB-ITEM FIELDS (A-ASSOCIATE-RQ and A-ASSOCIATE-AC)

### 27.252.2 Constructor & Destructor Documentation

27.252.2.1 `gdcm::network::SOPClassExtendedNegociationSub::SOPClassExtendedNegociationSub ( )`

### 27.252.3 Member Function Documentation

27.252.3.1 `void gdcm::network::SOPClassExtendedNegociationSub::Print ( std::ostream & os ) const`

27.252.3.2 `std::istream& gdcm::network::SOPClassExtendedNegociationSub::Read ( std::istream & is )`

27.252.3.3 `void gdcm::network::SOPClassExtendedNegociationSub::SetTuple ( const char * uid, uint8_t levelofsupport = 3, uint8_t levelofdigitalsig = 0, uint8_t elementcoercion = 2 )`

27.252.3.4 `size_t gdcm::network::SOPClassExtendedNegociationSub::Size ( ) const`

27.252.3.5 `const std::ostream& gdcm::network::SOPClassExtendedNegociationSub::Write ( std::ostream & os ) const`

The documentation for this class was generated from the following file:

- [gdcmSOPClassExtendedNegociationSub.h](#)

## 27.253 gdcm::SOPClassUIDToIOD Class Reference

Class convert a class SOP Class UID into [IOD](#).

```
#include <gdcmSOPClassUIDToIOD.h>
```

### Public Types

- typedef const char \* [const](#)(SOPClassUIDToIODType)[2]

### Static Public Member Functions

- static [const](#) char \* [GetIOD](#) (UIDs [const](#) &uid)
- static [const](#) char \* [GetIODFromSOPClassUID](#) ([const](#) char \*sopclassuid)
- static unsigned int [GetNumberOfSOPClassToIOD](#) ()  
*Return the number of SOP Class UID listed internally.*
- static [const](#) char \* [GetSOPClassUIDFromIOD](#) ([const](#) char \*iod)
- static SOPClassUIDToIODType & [GetSOPClassUIDToIOD](#) (unsigned int i)
- static SOPClassUIDToIODType \* [GetSOPClassUIDToIODs](#) ()

### 27.253.1 Detailed Description

Class convert a class SOP Class UID into [IOD](#).

Reference PS 3.4 [Table B.5-1 STANDARD SOP CLASSES](#)

### 27.253.2 Member Typedef Documentation

27.253.2.1 `typedef const char* gdcm::SOPClassUIDToIOD::const(SOPClassUIDToIODType)[2]`

### 27.253.3 Member Function Documentation

27.253.3.1 `static const char* gdcm::SOPClassUIDToIOD::GetIOD ( UIDs const & uid ) [static]`

Return the associated [IOD](#) based on a SOP Class UID uid (there is a one-to-one mapping from SOP Class UID to matching [IOD](#))

Examples:

[GenerateStandardSOPClasses.cxx](#).

27.253.3.2 `static const char* gdcm::SOPClassUIDToIOD::GetIODFromSOPClassUID ( const char * sopclassuid ) [static]`

27.253.3.3 `static unsigned int gdcm::SOPClassUIDToIOD::GetNumberOfSOPClassToIOD ( ) [static]`

Return the number of SOP Class UID listed internally.

27.253.3.4 `static const char* gdcm::SOPClassUIDToIOD::GetSOPClassUIDFromIOD ( const char * iod ) [static]`

27.253.3.5 `static SOPClassUIDToIODType& gdcm::SOPClassUIDToIOD::GetSOPClassUIDToIOD ( unsigned int i ) [static]`

27.253.3.6 `static SOPClassUIDToIODType* gdcm::SOPClassUIDToIOD::GetSOPClassUIDToIODs ( ) [static]`

The documentation for this class was generated from the following file:


- [gdcmSOPClassUIDToIOD.h](#)

## 27.254 gdcm::Sorter Class Reference


[Sorter](#) General class to do sorting using a custom function You simply need to provide a function of type: [Sorter::Sort←Function](#).

```
#include <gdcmSorter.h>
```

Inheritance diagram for gdcm::Sorter:

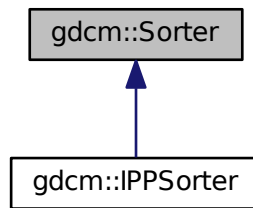

Collaboration diagram for gdcm::Sorter:

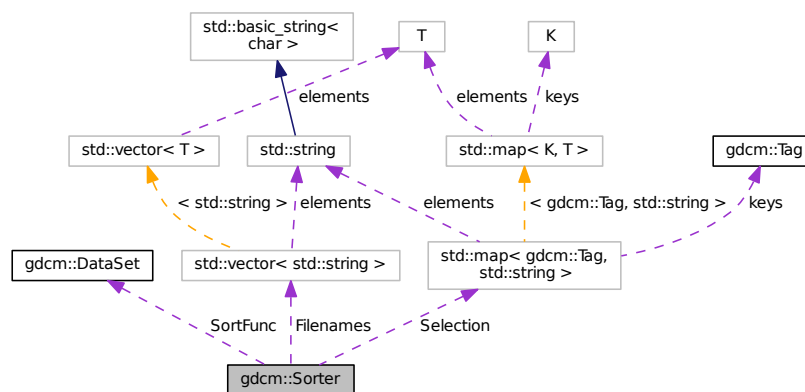

## Public Types

- typedef `bool(* SortFunction) (DataSet const &, DataSet const &)`  
Set the sort function which compares one dataset to the other.

## Public Member Functions

- `Sorter ()`
- virtual `~Sorter ()`
- `bool AddSelect (Tag const &tag, const char *value)`  
*UNSUPPORTED FOR NOW.*
- `const std::vector< std::string > & GetFilenames () const`
- `void Print (std::ostream &os) const`  
*Print.*

- void [SetSortFunction](#) ([SortFunction](#) f)
- virtual bool [Sort](#) (std::vector< std::string > const &filenames)  
*Typically the output of [Directory::GetFileNames\(\)](#)*
- virtual bool [StableSort](#) (std::vector< std::string > const &filenames)

## Protected Types

- typedef std::map< [Tag](#), std::string > [SelectionMap](#)

## Protected Attributes

- std::vector< std::string > [FileNames](#)
- std::map< [Tag](#), std::string > [Selection](#)
- [SortFunction](#) [SortFunc](#)

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Sorter](#) &s)

## 27.254.1 Detailed Description

[Sorter](#) General class to do sorting using a custom function You simply need to provide a function of type: [Sorter::SortFunction](#).

### Warning

implementation details. For now there is no cache mechanism. Which means that everytime you call Sort, all files specified as input paramater are *read*

### See also

[Scanner](#)

### Examples:

[SortImage.cxx](#), and [VolumeSorter.cxx](#).

## 27.254.2 Member Typedef Documentation

27.254.2.1 typedef std::map<[Tag](#),std::string> [gdcm::Sorter::SelectionMap](#) [protected]

27.254.2.2 typedef bool(\* [gdcm::Sorter::SortFunction](#)) ([DataSet](#) const &, [DataSet](#) const &)

Set the sort function which compares one dataset to the other.



### 27.254.3 Constructor & Destructor Documentation

27.254.3.1 `gdcm::Sorter::Sorter ( )`

27.254.3.2 `virtual gdcm::Sorter::~~Sorter ( )` `[virtual]`

### 27.254.4 Member Function Documentation

27.254.4.1 `bool gdcm::Sorter::AddSelect ( Tag const & tag, const char * value )`

UNSUPPORTED FOR NOW.

27.254.4.2 `const std::vector<std::string>& gdcm::Sorter::GetFileNames ( ) const` `[inline]`

Return the list of filenames as sorted by the specific algorithm used. Empty by default (before [Sort\(\)](#) is called)

Examples:

[Compute3DSpacing.cxx](#), [gdcmorthoplanes.cxx](#), [reslicesphere.cxx](#), [SortImage.cxx](#), and [VolumeSorter.cxx](#).

27.254.4.3 `void gdcm::Sorter::Print ( std::ostream & os ) const`

Print.

Examples:

[gdcmorthoplanes.cxx](#), [SortImage.cxx](#), and [VolumeSorter.cxx](#).

Referenced by `gdcm::operator<<()`.

27.254.4.4 `void gdcm::Sorter::SetSortFunction ( SortFunction f )`

Examples:

[SortImage.cxx](#), and [VolumeSorter.cxx](#).

27.254.4.5 `virtual bool gdcm::Sorter::Sort ( std::vector< std::string > const & filenames )` `[virtual]`

Typically the output of [Directory::GetFileNames\(\)](#)

Reimplemented in [gdcm::IPPSorter](#).

Examples:

[SortImage.cxx](#).

27.254.4.6 `virtual bool gdcm::Sorter::StableSort ( std::vector< std::string > const & filenames )` `[virtual]`

Examples:

[SortImage.cxx](#), and [VolumeSorter.cxx](#).

### 27.254.5 Friends And Related Function Documentation

27.254.5.1 `std::ostream& operator<< ( std::ostream & _os, const Sorter & s )` [*friend*]

### 27.254.6 Member Data Documentation

27.254.6.1 `std::vector<std::string> gdcm::Sorter::FileNames` [*protected*]

27.254.6.2 `std::map<Tag,std::string> gdcm::Sorter::Selection` [*protected*]

27.254.6.3 **SortFunction** `gdcm::Sorter::SortFunc` [*protected*]

The documentation for this class was generated from the following file:

- [gdcSorter.h](#)

## 27.255 **gdc**m::Spacing Class Reference

Class for [Spacing](#).

```
#include <gdcSpacing.h>
```

### Public Types

- enum [SpacingType](#) {  
[DETECTOR](#) = 0,  
[MAGNIFIED](#),  
[CALIBRATED](#),  
[UNKNOWN](#) }

### Public Member Functions

- [Spacing](#) ()
- [~Spacing](#) ()

### Static Public Member Functions

- static [Attribute](#)< 0x28, 0x34 > [ComputePixelAspectRatioFromPixelSpacing](#) (const [Attribute](#)< 0x28, 0x30 > &pixelspacing)

### 27.255.1 Detailed Description

Class for [Spacing](#).

It all began with a mail to WG6:

**Subject:** Imager Pixel [Spacing](#) vs Pixel [Spacing](#) **Body:** [Apologies for the duplicate post, namely to David Clunie & OFFIS team]

I have been trying to understand CP-586 in the following two cases:

On the one hand:

- DISCIMG/IMAGES/CRIMAGE taken from <http://dclunie.com/images/pixelspacingtestimages.zip> ↵

And on the other hand:

- [http://gdcm.sourceforge.net/thingies/cr\\_pixelspacing.dcm](http://gdcm.sourceforge.net/thingies/cr_pixelspacing.dcm)

If I understand correctly the CP, one is required to use Pixel [Spacing](#) for measurement ('true size' print) instead of Imager Pixel [Spacing](#), since the two attributes are present and Pixel [Spacing](#) is different from Imager Pixel [Spacing](#).

If this is correct, then the test data DISCIMG/IMAGES/CRIMAGE is incorrect. If this is incorrect (ie. I need to use Imager Pixel [Spacing](#)), then the display of cr\_pixelspacing.dcm for measurement will be incorrect.

Could someone please let me know what am I missing here? I could not find any information in any header that would allow me to differentiate those.

Thank you for your time,

Ref: <http://lists.nema.org/scripts/lyris.pl?sub=488573&id=400720477> See PS 3.3-2008, [Table C.7-11b IMAGE PIXEL MACRO ATTRIBUTES](#)

Ratio of the vertical size and horizontal size of the pixels in the image specified by a pair of integer values where the first value is the vertical pixel size, and the second value is the horizontal pixel size. Required if the aspect ratio values do not have a ratio of 1:1 and the physical pixel spacing is not specified by Pixel [Spacing](#) (0028,0030), or Imager Pixel [Spacing](#) (0018,1164) or Nominal Scanned Pixel [Spacing](#) (0018,2010), either for the entire [Image](#) or per-frame in a Functional Group [Macro](#). See C.7.6.3.1.7.

PS 3.3-2008 10.7.1.3 Pixel [Spacing Value](#) Order and Valid Values All pixel spacing related attributes shall have non-zero values, except when there is only a single row or column or pixel of data present, in which case the corresponding value may be zero.

Ref: [http://gdcm.sourceforge.net/wiki/index.php/Imager\\_Pixel\\_Spacing](http://gdcm.sourceforge.net/wiki/index.php/Imager_Pixel_Spacing)

## 27.255.2 Member Enumeration Documentation

### 27.255.2.1 enum gdcm::Spacing::SpacingType

Enumerator

**DETECTOR**  
**MAGNIFIED**  
**CALIBRATED**  
**UNKNOWN**

## 27.255.3 Constructor & Destructor Documentation

### 27.255.3.1 gdcm::Spacing::Spacing ( )

### 27.255.3.2 gdcm::Spacing::~~Spacing ( )

## 27.255.4 Member Function Documentation

### 27.255.4.1 static Attribute<0x28,0x34> gdcm::Spacing::ComputePixelAspectRatioFromPixelSpacing ( const Attribute< 0x28, 0x30 > & pixelspacing ) [static]

The documentation for this class was generated from the following file:

- [gdcmSpacing.h](#)

## 27.256 gdcm::Spectroscopy Class Reference

[Spectroscopy](#) class.

```
#include <gdcmSpectroscopy.h>
```

### Public Member Functions

- [Spectroscopy](#) ()

### 27.256.1 Detailed Description

[Spectroscopy](#) class.

### 27.256.2 Constructor & Destructor Documentation

27.256.2.1 [gdcm::Spectroscopy::Spectroscopy](#) ( ) [inline]

The documentation for this class was generated from the following file:

- [gdcmSpectroscopy.h](#)

## 27.257 gdcm::SplitMosaicFilter Class Reference

[SplitMosaicFilter](#) class Class to reshuffle bytes for a SIEMENS Mosaic image Siemens CSA [Image](#) Header CSA:= Common Siemens Architecture, sometimes also known as Common syngo Architecture.

```
#include <gdcmSplitMosaicFilter.h>
```

### Public Member Functions

- [SplitMosaicFilter](#) ()
- [~SplitMosaicFilter](#) ()
- bool [ComputeMOSAICDimensions](#) (unsigned int dims[3])
- [File](#) & [GetFile](#) ()
- const [File](#) & [GetFile](#) () const
- const [Image](#) & [GetImage](#) () const
- [Image](#) & [GetImage](#) ()
- void [SetFile](#) (const [File](#) &f)
- void [SetImage](#) (const [Image](#) &image)
- bool [Split](#) ()

*Split the SIEMENS MOSAIC image.*

### 27.257.1 Detailed Description

[SplitMosaicFilter](#) class Class to reshuffle bytes for a SIEMENS Mosaic image Siemens CSA [Image](#) Header CSA:= Common Siemens Architecture, sometimes also known as Common syngo Architecture.

### 27.257.2 Constructor & Destructor Documentation

27.257.2.1 `gdcm::SplitMosaicFilter::SplitMosaicFilter ( )`

27.257.2.2 `gdcm::SplitMosaicFilter::~~SplitMosaicFilter ( )`

### 27.257.3 Member Function Documentation

27.257.3.1 `bool gdcm::SplitMosaicFilter::ComputeMOSAICDimensions ( unsigned int dims[3] )`

Compute the new dimensions according to private information stored in the MOSAIC header.

27.257.3.2 `File& gdcm::SplitMosaicFilter::GetFile ( )` `[inline]`

27.257.3.3 `const File& gdcm::SplitMosaicFilter::GetFile ( ) const` `[inline]`

27.257.3.4 `const Image& gdcm::SplitMosaicFilter::GetImage ( ) const` `[inline]`

27.257.3.5 `Image& gdcm::SplitMosaicFilter::GetImage ( )` `[inline]`

27.257.3.6 `void gdcm::SplitMosaicFilter::SetFile ( const File & f )` `[inline]`

27.257.3.7 `void gdcm::SplitMosaicFilter::SetImage ( const Image & image )`

27.257.3.8 `bool gdcm::SplitMosaicFilter::Split ( )`

Split the SIEMENS MOSAIC image.


The documentation for this class was generated from the following file:

- [gdcmSplitMosaicFilter.h](#)


## 27.258 gdcm::StartEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::StartEvent`:

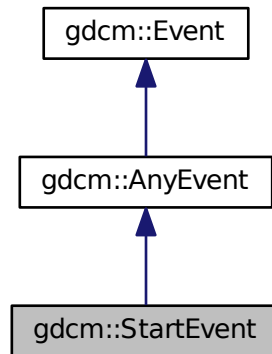

Collaboration diagram for `gdcm::StartEvent`:

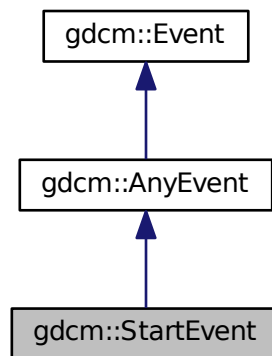

### Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmEvent.h](#)

## 27.259 `gdcm::static_assert_test< x >` Struct Template Reference

```
#include <gdcmStaticAssert.h>
```

The documentation for this struct was generated from the following file:

- [gdcmStaticAssert.h](#)

## 27.260 `gdcm::STATIC_ASSERTION_FAILURE< x >` Struct Template Reference

```
#include <gdcmStaticAssert.h>
```

The documentation for this struct was generated from the following file:

- [gdcmStaticAssert.h](#)

## 27.261 `gdcm::STATIC_ASSERTION_FAILURE< true >` Struct Template Reference

```
#include <gdcmStaticAssert.h>
```

### Public Types

- enum { [value](#) = 1 }

### 27.261.1 Member Enumeration Documentation

#### 27.261.1.1 anonymous enum

Enumerator

***value***

The documentation for this struct was generated from the following file:

- [gdcmStaticAssert.h](#)

## 27.262 `gdcm::StreamImageReader` Class Reference

[StreamImageReader](#).

```
#include <gdcmStreamImageReader.h>
```

### Public Member Functions

- [StreamImageReader](#) ()
- virtual [~StreamImageReader](#) ()
- bool [CanReadImage](#) () const

- void [DefinePixelExtent](#) (uint16\_t inXMin, uint16\_t inXMax, uint16\_t inYMin, uint16\_t inYMax, uint16\_t inZMin=0, uint16\_t inZMax=1)
- uint32\_t [DefineProperBufferLength](#) () const
- std::vector< unsigned int > [GetDimensionsValueForResolution](#) (unsigned int)
- [File](#) const & [GetFile](#) () const
- bool [Read](#) (char \*inReadBuffer, const std::size\_t &inBufferLength)
- virtual bool [ReadImageInformation](#) ()
- void [SetFileName](#) (const char \*inFileName)
- void [SetStream](#) (std::istream &inStream)

### 27.262.1 Detailed Description

[StreamImageReader](#).

#### Note

its role is to convert the DICOM [DataSet](#) into a [Image](#) representation via an ITK streaming (ie, multithreaded) interface [Image](#) is different from [Pixmap](#) has it has a position and a direction in Space. Currently, this class is thread safe in that it can read a single extent in a single thread. Multiple versions can be used for multiple extents/threads.

See also

[Image](#)

Examples:

[StreamImageReaderTest.cxx](#).

### 27.262.2 Constructor & Destructor Documentation

27.262.2.1 `gdcm::StreamImageReader::StreamImageReader ( )`

27.262.2.2 `virtual gdcm::StreamImageReader::~~StreamImageReader ( ) [virtual]`

### 27.262.3 Member Function Documentation

27.262.3.1 `bool gdcm::StreamImageReader::CanReadImage ( ) const`

Only RAW images are currently readable by the stream reader. As more streaming codecs are added, then this function will be updated to reflect those changes. Calling this function prior to reading will ensure that only streamable files are streamed. Make sure to call [ReadImageInformation](#) prior to calling this function.

Examples:

[StreamImageReaderTest.cxx](#).

27.262.3.2 `void gdcm::StreamImageReader::DefinePixelExtent ( uint16_t inXMin, uint16_t inXMax, uint16_t inYMin, uint16_t inYMax, uint16_t inZMin = 0, uint16_t inZMax = 1 )`

Defines an image extent for the [Read](#) function. DICOM states that an image can have no more than  $2^{16}$  pixels per edge (as of 2009) In this case, the pixel extents ignore the direction cosines entirely, and assumes that the origin of the



image is at location 0,0 (regardless of the definition in space per the tags). So, if the first 100 pixels of the first row are to be read in, this function should be called with DefinePixelExtent(0, 100, 0, 1), regardless of pixel size or orientation.

Examples:

[StreamImageReaderTest.cxx](#).

**27.262.3.3** `uint32_t gdcm::StreamImageReader::DefineProperBufferLength ( ) const`

Paying attention to the pixel format and so forth, define the proper buffer length for the user. The return amount is in bytes. Call this function to determine the size of the char\* buffer that will need to be passed in to ReadImageSubregion(). If the return is 0, then that means that the pixel extent was not defined prior

Examples:

[StreamImageReaderTest.cxx](#).

**27.262.3.4** `std::vector<unsigned int> gdcm::StreamImageReader::GetDimensionsValueForResolution ( unsigned int )`

**27.262.3.5** `File const& gdcm::StreamImageReader::GetFile ( ) const`

Returns the dataset read by ReadImageInformation Couple this with the [ImageHelper](#) to get statistics about the image, like pixel extent, to be able to initialize buffers for reading

Examples:

[StreamImageReaderTest.cxx](#).

**27.262.3.6** `bool gdcm::StreamImageReader::Read ( char * inReadBuffer, const std::size_t & inBufferLength )`

Read the DICOM image. There are three reasons for failure:

1. The extent is not set
2. the conversion from char\* to std::ostream (internally) fails
3. the given buffer isn't large enough to accommodate the desired pixel extent. This method has been implemented to look similar to the metainageio in itk MUST have an extent defined, or else Read will return false. If no particular extent is required, use [ImageReader](#) instead.

Examples:

[StreamImageReaderTest.cxx](#).

**27.262.3.7** `virtual bool gdcm::StreamImageReader::ReadImageInformation ( ) [virtual]`

Set the spacing and dimension information for the set filename. returns false if the file is not initialized or not an image, with the pixel (7fe0,0010) tag.

Examples:

[StreamImageReaderTest.cxx](#).

27.262.3.8 void gdcM::StreamImageReader::SetFileName ( const char \* *inFileName* )

One of either SetFileName or SetStream must be called prior to any other functions. These initialize an internal [Reader](#) class to be able to get non-pixel image information.

Examples:

[ExtractOneFrame.cs](#), and [StreamImageReaderTest.cxx](#).

27.262.3.9 void gdcM::StreamImageReader::SetStream ( std::istream & *inStream* )

The documentation for this class was generated from the following file:


- [gdcMStreamImageReader.h](#)

## 27.263 gdcM::StreamImageWriter Class Reference

[StreamImageReader](#).

```
#include <gdcMStreamImageWriter.h>
```

Collaboration diagram for gdcm::StreamImageWriter:

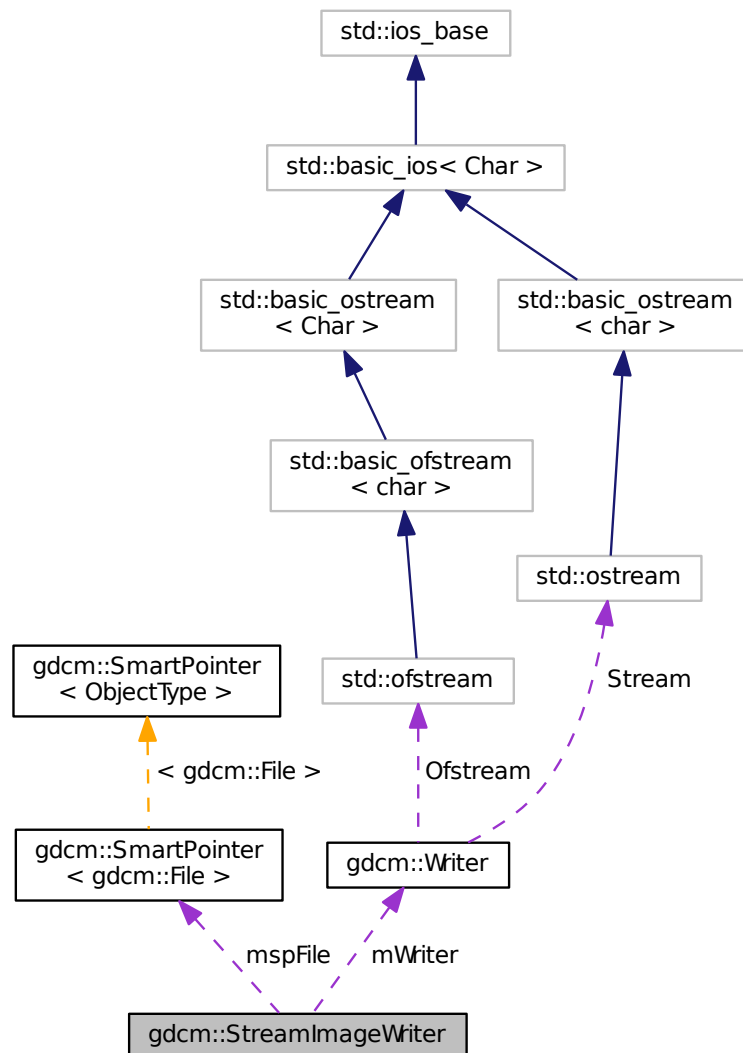

## Public Member Functions

- [StreamImageWriter](#) ()
- virtual [~StreamImageWriter](#) ()
- bool [CanWriteFile](#) () const
- void [DefinePixelExtent](#) (uint16\_t inXMin, uint16\_t inXMax, uint16\_t inYMin, uint16\_t inYMax, uint16\_t inZMin=0, uint16\_t inZMax=1)
- uint32\_t [DefineProperBufferLength](#) ()
- void [SetFile](#) (const [File](#) &inFile)
- void [SetFileName](#) (const char \*inFileName)

- void [SetStream](#) (std::ostream &inStream)
- bool [Write](#) (void \*inWriteBuffer, const std::size\_t &inBufferLength)
- virtual bool [WriteImageInformation](#) ()

### Protected Member Functions

- virtual bool [WriteImageSubregionRAW](#) (char \*inWriteBuffer, const std::size\_t &inBufferLength)
- int [WriteRawHeader](#) ([RAWCodec](#) \*inCodec, std::ostream \*inStream)

### Protected Attributes

- int [mElementOffsets](#)
- int [mElementOffsets1](#)
- [SmartPointer](#)< [File](#) > [mspFile](#)
- [Writer](#) [mWriter](#)
- uint16\_t [mXMax](#)
- uint16\_t [mXMin](#)
- uint16\_t [mYMax](#)
- uint16\_t [mYMin](#)
- uint16\_t [mZMax](#)
- uint16\_t [mZMin](#)

## 27.263.1 Detailed Description

[StreamImageReader](#).

### Note

its role is to convert the DICOM [DataSet](#) into a [Image](#) representation via an ITK streaming (ie, multithreaded) interface [Image](#) is different from [Pixmap](#) has it has a position and a direction in Space. Currently, this class is threadsafe in that it can read a single extent in a single thread. Multiple versions can be used for multiple extents/threads.

### See also

[Image](#)

### Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

## 27.263.2 Constructor & Destructor Documentation

27.263.2.1 `gdcm::StreamImageWriter::StreamImageWriter ( )`

27.263.2.2 `virtual gdcm::StreamImageWriter::~~StreamImageWriter ( )` `[virtual]`

## 27.263.3 Member Function Documentation

### 27.263.3.1 bool gdcm::StreamImageWriter::CanWriteFile ( ) const

This function determines if a file can even be written using the streaming writer unlike the reader, can be called before WriteImageInformation, but must be called after SetFile.

Examples:

[Extracting\\_All\\_Resolution.cxx](#), and [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#).

### 27.263.3.2 void gdcm::StreamImageWriter::DefinePixelExtent ( uint16\_t *inXMin*, uint16\_t *inXMax*, uint16\_t *inYMin*, uint16\_t *inYMax*, uint16\_t *inZMin* = 0, uint16\_t *inZMax* = 1 )

Defines an image extent for the Read function. DICOM states that an image can have no more than  $2^{16}$  pixels per edge (as of 2009) In this case, the pixel extents ignore the direction cosines entirely, and assumes that the origin of the image is at location 0,0 (regardless of the definition in space per the tags). So, if the first 100 pixels of the first row are to be read in, this function should be called with DefinePixelExtent(0, 100, 0, 1), regardless of pixel size or orientation.  
15 nov 2010: added z dimension, defaults to being 1 plane large

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

### 27.263.3.3 uint32\_t gdcm::StreamImageWriter::DefineProperBufferLength ( )

Paying attention to the pixel format and so forth, define the proper buffer length for the user. The return amount is in bytes. If the return is 0, then that means that the pixel extent was not defined prior this return is for RAW inputs which are then encoded by the writer, but are used to ensure that the writer gets the proper buffer size

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

### 27.263.3.4 void gdcm::StreamImageWriter::SetFile ( const File & *inFile* )

Set the image information to be written to disk that is everything but the pixel information: (7fe0,0010) PixelData

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

### 27.263.3.5 void gdcm::StreamImageWriter::SetFileName ( const char \* *inFileName* )

One of either SetFileName or SetStream must be called prior to any other functions. These initialize an internal [Reader](#) class to be able to get non-pixel image information.

### 27.263.3.6 void gdcm::StreamImageWriter::SetStream ( std::ostream & *inStream* )

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

**27.263.3.7** `bool gdcmm::StreamImageWriter::Write ( void * inWriteBuffer, const std::size_t & inBufferLength )`

Read the DICOM image. There are three reasons for failure:

1. The extent is not set
2. the conversion from void\* to std::ostream (internally) fails
3. the given buffer isn't large enough to accomodate the desired pixel extent. This method has been implemented to look similar to the metaimageio in itk MUST have an extent defined, or else Read will return false. If no particular extent is required, use [ImageReader](#) instead.

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

**27.263.3.8** `virtual bool gdcmm::StreamImageWriter::WriteImageInformation ( ) [virtual]`

Write the header information to disk, and a bunch of zeros for the actual pixel information. Of course, if we're doing a non-compressed format, that works but if it's compressed, we have to force the ordering of chunks that are written.

Examples:

[Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

**27.263.3.9** `virtual bool gdcmm::StreamImageWriter::WriteImageSubregionRAW ( char * inWriteBuffer, const std::size_t & inBufferLength ) [protected], [virtual]`

Using the min, max, etc set by DefinePixelExtent, this will fill the given buffer. Make sure to call DefinePixelExtent and to initialize the buffer with the amount given by DefineProperBufferLength prior to calling this. reads by the RAW codec; other codecs are added once implemented

**27.263.3.10** `int gdcmm::StreamImageWriter::WriteRawHeader ( RAWCodec * inCodec, std::ostream * inStream ) [protected]`

when writing a raw file, we know the full extent, and can just write the first 12 bytes out (the tag, the [VR](#), and the size) when we do compressed files, we'll do it in chunks, as described in 2009-3, part 5, Annex A, section 4. Pass the raw codec so that in the rare case of a bigendian explicit raw, the first 12 bytes written out should still be kosher. returns -1 if there's any failure, or the complete offset (12 bytes) if it works. Those 12 bytes are then added to the position in order to determine where to write.

## 27.263.4 Member Data Documentation

**27.263.4.1** `int gdcmm::StreamImageWriter::mElementOffsets [protected]`

The result of WriteRawHeader (or another header, when that's implemented) This result is saved so that the first N bytes aren't constantly being rewritten for each chunk that's passed in. For compressed data, the offset table will require rewrites of data.

27.263.4.2 int gdcm::StreamImageWriter::mElementOffsets1 [protected]

27.263.4.3 SmartPointer<File> gdcm::StreamImageWriter::mspFile [protected]

27.263.4.4 Writer gdcm::StreamImageWriter::mWriter [protected]

27.263.4.5 uint16\_t gdcm::StreamImageWriter::mXMax [protected]

27.263.4.6 uint16\_t gdcm::StreamImageWriter::mXMin [protected]

27.263.4.7 uint16\_t gdcm::StreamImageWriter::mYMax [protected]

27.263.4.8 uint16\_t gdcm::StreamImageWriter::mYMin [protected]

27.263.4.9 uint16\_t gdcm::StreamImageWriter::mZMax [protected]

27.263.4.10 uint16\_t gdcm::StreamImageWriter::mZMin [protected]

The documentation for this class was generated from the following file:


- [gdcmStreamImageWriter.h](#)

## 27.264 gdcm::String< TDelimiter, TMaxLength, TPadChar > Class Template Reference


[String.](#)

```
#include <gdcmString.h>
```

Inheritance diagram for gdcm::String< TDelimiter, TMaxLength, TPadChar >:

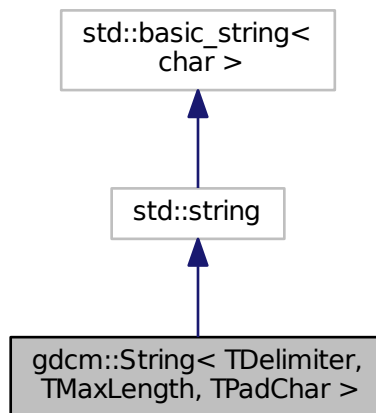

Collaboration diagram for `gdcm::String< TDelimiter, TMaxLength, TPadChar >`:

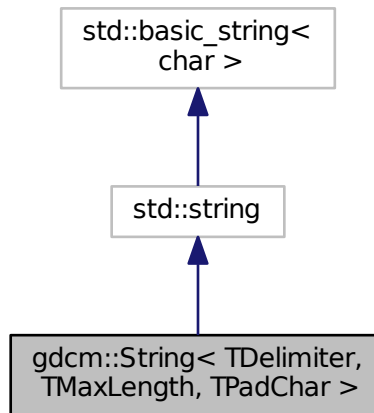

## Public Types

- typedef `std::string::const_iterator` [const\\_iterator](#)
- typedef `std::string::const_reference` [const\\_reference](#)
- typedef `std::string::const_reverse_iterator` [const\\_reverse\\_iterator](#)
- typedef `std::string::difference_type` [difference\\_type](#)
- typedef `std::string::iterator` [iterator](#)
- typedef `std::string::pointer` [pointer](#)
- typedef `std::string::reference` [reference](#)
- typedef `std::string::reverse_iterator` [reverse\\_iterator](#)
- typedef `std::string::size_type` [size\\_type](#)
- typedef `std::string::value_type` [value\\_type](#)

## Public Member Functions

- [String](#) ()  
*String constructors.*
- [String](#) (const [value\\_type](#) \*s)
- [String](#) (const [value\\_type](#) \*s, [size\\_type](#) n)
- [String](#) (const std::string &s, [size\\_type](#) pos=0, [size\\_type](#) n=npos)
- bool [IsValid](#) () const  
*return if string is valid*
- [operator const char \\*](#) () const  
*WARNING: Trailing \0 might be lost in this operation:*
- std::string [Trim](#) () const
- [gdcm::String](#)< TDelimiter, TMaxLength, TPadChar > [Truncate](#) () const



## Static Public Member Functions

- static std::string [Trim](#) (const char \*input)

### 27.264.1 Detailed Description

```
template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '>class gdcmm::String< TDelimiter, TMaxLength, TPadChar >
```

[String](#).

#### Note

TDelimiter template parameter is used to separate multiple [String](#) (VM1 >) TMaxLength is only a hint. Noone actually respect the max length TPadChar is the string padding (0 or space)

### 27.264.2 Member Typedef Documentation

27.264.2.1 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::const_iterator gdcmm::String< TDelimiter, TMaxLength, TPadChar >::const_iterator`

27.264.2.2 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::const_reference gdcmm::String< TDelimiter, TMaxLength, TPadChar >::const_reference`

27.264.2.3 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::const_reverse_iterator gdcmm::String< TDelimiter, TMaxLength, TPadChar >::const_reverse_iterator`

27.264.2.4 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::difference_type gdcmm::String< TDelimiter, TMaxLength, TPadChar >::difference_type`

27.264.2.5 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::iterator gdcmm::String< TDelimiter, TMaxLength, TPadChar >::iterator`

27.264.2.6 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::pointer gdcmm::String< TDelimiter, TMaxLength, TPadChar >::pointer`

27.264.2.7 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::reference gdcmm::String< TDelimiter, TMaxLength, TPadChar >::reference`

27.264.2.8 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::reverse_iterator gdcmm::String< TDelimiter, TMaxLength, TPadChar >::reverse_iterator`

27.264.2.9 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::size_type gdcmm::String< TDelimiter, TMaxLength, TPadChar >::size_type`

27.264.2.10 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = ' '> typedef std::string::value_type gdcmm::String< TDelimiter, TMaxLength, TPadChar >::value_type`

### 27.264.3 Constructor & Destructor Documentation

27.264.3.1 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar >::String ( ) [inline]`

[String](#) constructors.

27.264.3.2 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar >::String ( const value_type * s ) [inline]`

27.264.3.3 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar >::String ( const value_type * s, size_type n ) [inline]`

27.264.3.4 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar >::String ( const std::string & s, size_type pos = 0, size_type n = npos ) [inline]`

## 27.264.4 Member Function Documentation

27.264.4.1 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> bool gdcmm::String< TDelimiter, TMaxLength, TPadChar >::IsValid ( ) const [inline]`

return if string is valid

Referenced by `gdcmm::String< TDelimiter, TMaxLength, TPadChar >::Truncate()`.

27.264.4.2 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar >::operator const char * ( ) const [inline]`

WARNING: Trailing \0 might be lost in this operation:

27.264.4.3 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> std::string gdcmm::String< TDelimiter, TMaxLength, TPadChar >::Trim ( ) const [inline]`

Trim function is required to return a `std::string` object, otherwise we could not create a [gdcmm::String](#) object with an odd number of bytes...

27.264.4.4 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> static std::string gdcmm::String< TDelimiter, TMaxLength, TPadChar >::Trim ( const char * input ) [inline], [static]`

27.264.4.5 `template<char TDelimiter = '\\', unsigned int TMaxLength = 64, char TPadChar = '> gdcmm::String< TDelimiter, TMaxLength, TPadChar > gdcmm::String< TDelimiter, TMaxLength, TPadChar >::Truncate ( ) const [inline]`

References `gdcmm::String< TDelimiter, TMaxLength, TPadChar >::IsValid()`.

The documentation for this class was generated from the following file:

- [gdcmmString.h](#)

## 27.265 gdcmm::StringFilter Class Reference

[StringFilter](#) [StringFilter](#) is the class that make `gdcmm2.x` looks more like `gdcmm1` and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.

```
#include <gdcmStringFilter.h>
```

## Public Member Functions

- [StringFilter](#) ()
- [~StringFilter](#) ()
- bool [ExecuteQuery](#) (std::string const &query, std::string &value) const
- std::string [FromString](#) (const [Tag](#) &t, const char \*value, [VL](#) const &vl)
- std::string [FromString](#) (const [Tag](#) &t, const char \*value, size\_t len)  
*Convert to string the char array defined by the pair (value,len)*
- [File](#) & [GetFile](#) ()
- const [File](#) & [GetFile](#) () const
- void [SetDicts](#) (const [Dicts](#) &dicts)  
*Allow user to pass in there own dicts.*
- void [SetFile](#) (const [File](#) &f)  
*Set/Get File.*
- std::string [ToString](#) (const [DataElement](#) &de) const
- std::string [ToString](#) (const [Tag](#) &t) const  
*Directly from a Tag:*
- std::pair< std::string, std::string > [ToStringPair](#) (const [DataElement](#) &de) const
- std::pair< std::string, std::string > [ToStringPair](#) (const [Tag](#) &t) const  
*Directly from a Tag:*
- void [UseDictAlways](#) (bool)

## Protected Member Functions

- bool [ExecuteQuery](#) (std::string const &query, [DataSet](#) const &ds, std::string &value) const
- std::pair< std::string, std::string > [ToStringPair](#) (const [Tag](#) &t, [DataSet](#) const &ds) const

### 27.265.1 Detailed Description

[StringFilter](#) [StringFilter](#) is the class that make gdcm2.x looks more like gdcm1 and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.

Examples:

[ReadAndPrintAttributes.cxx](#).

### 27.265.2 Constructor & Destructor Documentation

27.265.2.1 [gdcm::StringFilter::StringFilter](#) ( )

27.265.2.2 [gdcm::StringFilter::~~StringFilter](#) ( )

### 27.265.3 Member Function Documentation

27.265.3.1 bool [gdcm::StringFilter::ExecuteQuery](#) ( std::string const & *query*, std::string & *value* ) const

Execute the XPATH query to find a value (as string) return false when attribute is not found (or an error in the XPATH query) You need to make sure that your XPATH query is syntatically correct

27.265.3.2 `bool gdcM::StringFilter::ExecuteQuery ( std::string const & query, DataSet const & ds, std::string & value ) const`  
`[protected]`

27.265.3.3 `std::string gdcM::StringFilter::FromString ( const Tag & t, const char * value, VL const & vl )`

27.265.3.4 `std::string gdcM::StringFilter::FromString ( const Tag & t, const char * value, size_t len )`

Convert to string the char array defined by the pair (value,len)

27.265.3.5 `File& gdcM::StringFilter::GetFile ( ) [inline]`

27.265.3.6 `const File& gdcM::StringFilter::GetFile ( ) const [inline]`

27.265.3.7 `void gdcM::StringFilter::SetDicts ( const Dicts & dicts )`

Allow user to pass in there own dicts.

27.265.3.8 `void gdcM::StringFilter::SetFile ( const File & f ) [inline]`

Set/Get [File](#).

Examples:

[ReadAndPrintAttributes.cxx](#), and [SimplePrintPatientName.cs](#).

27.265.3.9 `std::string gdcM::StringFilter::ToString ( const DataElement & de ) const`

Convert to string the [ByteValue](#) contained in a [DataElement](#). The [DataElement](#) must be coming from the actual [DataSet](#) associated with [File](#) (see [SetFile](#)).

Examples:

[ReadAndPrintAttributes.cxx](#).

27.265.3.10 `std::string gdcM::StringFilter::ToString ( const Tag & t ) const`

Directly from a [Tag](#):

27.265.3.11 `std::pair<std::string, std::string> gdcM::StringFilter::ToStringPair ( const DataElement & de ) const`

Convert to string the [ByteValue](#) contained in a [DataElement](#) the returned elements are: pair.first : the name as found in the dictionary of [DataElement](#) pair.second : the value encoded into a string (US,UL...) are properly converted

Examples:

[ReadAndPrintAttributes.cxx](#).

27.265.3.12 `std::pair<std::string, std::string> gdcM::StringFilter::ToStringPair ( const Tag & t ) const`

Directly from a [Tag](#):

27.265.3.13 `std::pair<std::string, std::string> gdcm::StringFilter::ToStringPair ( const Tag & t, DataSet const & ds ) const`  
[protected]

27.265.3.14 `void gdcm::StringFilter::UseDictAlways ( bool )` [inline]

The documentation for this class was generated from the following file:

- [gdcmStringFilter.h](#)

## 27.266 gdcm::Study Class Reference

[Study.](#)

```
#include <gdcmStudy.h>
```

### Public Member Functions

- [Study \(\)](#)

### 27.266.1 Detailed Description

[Study.](#)

### 27.266.2 Constructor & Destructor Documentation

27.266.2.1 `gdcm::Study::Study ( )` [inline]

The documentation for this class was generated from the following file:


- [gdcmStudy.h](#)

## 27.267 gdcm::Subject Class Reference


[Subject.](#)

```
#include <gdcmSubject.h>
```

Inheritance diagram for `gdcmm::Subject`:

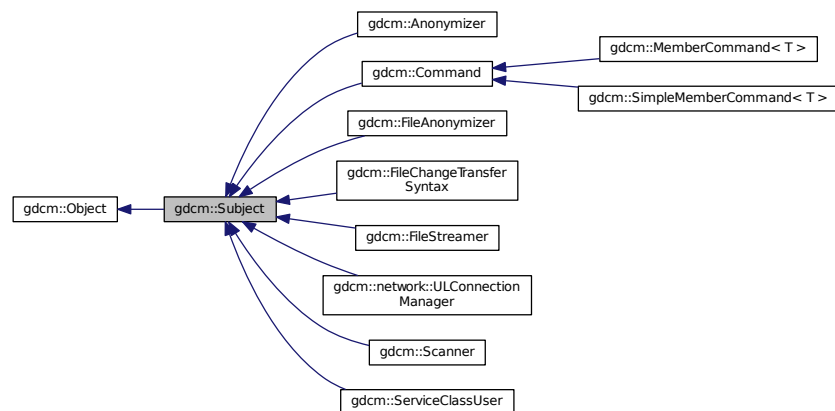

Collaboration diagram for `gdcmm::Subject`:

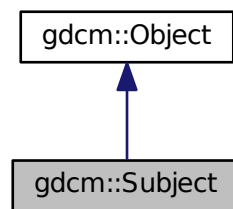

## Public Member Functions

- [Subject](#) ()
- [~Subject](#) ()
- unsigned long [AddObserver](#) (const [Event](#) &event, [Command](#) \*)
- unsigned long [AddObserver](#) (const [Event](#) &event, [Command](#) \*) const
- [Command](#) \* [GetCommand](#) (unsigned long tag)
- bool [HasObserver](#) (const [Event](#) &event) const
- void [InvokeEvent](#) (const [Event](#) &)
- void [InvokeEvent](#) (const [Event](#) &) const
- void [RemoveAllObservers](#) ()
- void [RemoveObserver](#) (unsigned long tag)

## Additional Inherited Members

### 27.267.1 Detailed Description

[Subject](#).

See also

[Command Event](#)

Examples:

[SimpleScanner.cxx](#).

### 27.267.2 Constructor & Destructor Documentation

27.267.2.1 `gdcmm::Subject::Subject ( )`

27.267.2.2 `gdcmm::Subject::~~Subject ( )`

### 27.267.3 Member Function Documentation

27.267.3.1 `unsigned long gdcmm::Subject::AddObserver ( const Event & event, Command * )`

Allow people to add/remove/invoke observers (callbacks) to any GDCM object. This is an implementation of the subject/observer design pattern. An observer is added by specifying an event to respond to and an [gdcmm::Command](#) to execute. It returns an unsigned long tag which can be used later to remove the event or retrieve the command. The memory for the [Command](#) becomes the responsibility of this object, so don't pass the same instance of a command to two different objects

27.267.3.2 `unsigned long gdcmm::Subject::AddObserver ( const Event & event, Command * ) const`

27.267.3.3 `Command* gdcmm::Subject::GetCommand ( unsigned long tag )`

Get the command associated with the given tag. NOTE: This returns a pointer to a [Command](#), but it is safe to assign this to a `Command::Pointer`. Since [Command](#) inherits from `LightObject`, at this point in the code, only a pointer or a reference to the [Command](#) can be used.

27.267.3.4 `bool gdcmm::Subject::HasObserver ( const Event & event ) const`

Return true if an observer is registered for this event.

27.267.3.5 `void gdcmm::Subject::InvokeEvent ( const Event & )`

Call `Execute` on all the `Commands` observing this event id.

27.267.3.6 `void gdcmm::Subject::InvokeEvent ( const Event & ) const`

Call `Execute` on all the `Commands` observing this event id. The actions triggered by this call doesn't modify this object.

27.267.3.7 void `gdcm::Subject::RemoveAllObservers` ( )

Remove all observers .

27.267.3.8 void `gdcm::Subject::RemoveObserver` ( unsigned long *tag* )

Remove the observer with this tag value.

The documentation for this class was generated from the following file:


- [gdcmSubject.h](#)

## 27.268 `gdcm::Surface` Class Reference


This class defines a SURFACE IE. This members are taken from required surface mesh module attributes.

```
#include <gdcmSurface.h>
```

Inheritance diagram for `gdcm::Surface`:

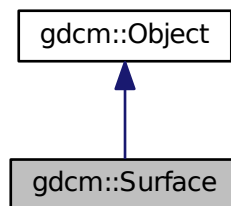

Collaboration diagram for `gdcm::Surface`:

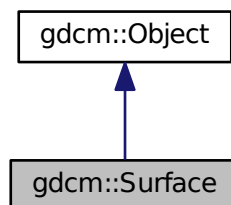



## Public Types

- enum STATES {  
NO = 0,  
YES,  
UNKNOWN,  
STATES\_END }
- enum VIEWType {  
SURFACE = 0,  
WIREFRAME,  
POINTS,  
VIEWType\_END }

Enumeration for Recommended Presentation [Type](#).

## Public Member Functions

- [Surface](#) ()
- virtual [~Surface](#) ()
- [SegmentHelper::BasicCodedEntry](#) const & [GetAlgorithmFamily](#) () const
- [SegmentHelper::BasicCodedEntry](#) & [GetAlgorithmFamily](#) ()
- const char \* [GetAlgorithmName](#) () const
- const char \* [GetAlgorithmVersion](#) () const
- const float \* [GetAxisOfRotation](#) () const
- const float \* [GetCenterOfRotation](#) () const
- STATES [GetFiniteVolume](#) () const
- STATES [GetManifold](#) () const
- float [GetMaximumPointDistance](#) () const
- float [GetMeanPointDistance](#) () const
- [MeshPrimitive](#) const & [GetMeshPrimitive](#) () const
- [MeshPrimitive](#) & [GetMeshPrimitive](#) ()
- unsigned long [GetNumberOfSurfacePoints](#) () const
- unsigned long [GetNumberOfVectors](#) () const
- const [DataElement](#) & [GetPointCoordinatesData](#) () const
- [DataElement](#) & [GetPointCoordinatesData](#) ()
- const float \* [GetPointPositionAccuracy](#) () const
- const float \* [GetPointsBoundingBoxCoordinates](#) () const
- [SegmentHelper::BasicCodedEntry](#) const & [GetProcessingAlgorithm](#) () const
- [SegmentHelper::BasicCodedEntry](#) & [GetProcessingAlgorithm](#) ()
- const unsigned short \* [GetRecommendedDisplayCIELabValue](#) () const
- unsigned short [GetRecommendedDisplayCIELabValue](#) (const unsigned int idx) const
- unsigned short [GetRecommendedDisplayGrayscaleValue](#) () const
- float [GetRecommendedPresentationOpacity](#) () const
- VIEWType [GetRecommendedPresentationType](#) () const
- const char \* [GetSurfaceComments](#) () const
- unsigned long [GetSurfaceNumber](#) () const
- bool [GetSurfaceProcessing](#) () const
- const char \* [GetSurfaceProcessingDescription](#) () const
- float [GetSurfaceProcessingRatio](#) () const
- const float \* [GetVectorAccuracy](#) () const
- const [DataElement](#) & [GetVectorCoordinateData](#) () const
- [DataElement](#) & [GetVectorCoordinateData](#) ()

- unsigned short [GetVectorDimensionality](#) () const
- void [SetAlgorithmFamily](#) ([SegmentHelper::BasicCodedEntry](#) const &BSE)
- void [SetAlgorithmName](#) (const char \*str)
- void [SetAlgorithmVersion](#) (const char \*str)
- void [SetAxisOfRotation](#) (const float \*axis)
- void [SetCenterOfRotation](#) (const float \*center)
- void [SetFiniteVolume](#) ([STATES](#) state)
- void [SetManifold](#) ([STATES](#) state)
- void [SetMaximumPointDistance](#) (float maximum)
- void [SetMeanPointDistance](#) (float average)
- void [SetMeshPrimitive](#) ([MeshPrimitive](#) &mp)
- void [SetNumberOfSurfacePoints](#) (const unsigned long nb)
- void [SetNumberOfVectors](#) (const unsigned long nb)
- void [SetPointCoordinatesData](#) ([DataElement](#) const &de)
- void [SetPointPositionAccuracy](#) (const float \*accuracies)
- void [SetPointsBoundingBoxCoordinates](#) (const float \*coordinates)
- void [SetProcessingAlgorithm](#) ([SegmentHelper::BasicCodedEntry](#) const &BSE)
- void [SetRecommendedDisplayCIELabValue](#) (const unsigned short vl[3])
- void [SetRecommendedDisplayCIELabValue](#) (const unsigned short vl, const unsigned int idx=0)
- void [SetRecommendedDisplayCIELabValue](#) (const std::vector< unsigned short > &vl)
- void [SetRecommendedDisplayGrayscaleValue](#) (const unsigned short vl)
- void [SetRecommendedPresentationOpacity](#) (const float opacity)
- void [SetRecommendedPresentationType](#) ([VIEWType](#) type)
- void [SetSurfaceComments](#) (const char \*comment)
- void [SetSurfaceNumber](#) (const unsigned long nb)
- void [SetSurfaceProcessing](#) (bool b)
- void [SetSurfaceProcessingDescription](#) (const char \*description)
- void [SetSurfaceProcessingRatio](#) (const float ratio)
- void [SetVectorAccuracy](#) (const float \*accuracy)
- void [SetVectorCoordinateData](#) ([DataElement](#) const &de)
- void [SetVectorDimensionality](#) (const unsigned short dim)

## Static Public Member Functions

- static [STATES](#) [GetSTATES](#) (const char \*state)
- static const char \* [GetSTATESString](#) ([STATES](#) state)
- static [VIEWType](#) [GetVIEWType](#) (const char \*type)
- static const char \* [GetVIEWTypeString](#) ([VIEWType](#) type)

## Additional Inherited Members

### 27.268.1 Detailed Description

This class defines a SURFACE IE. This members are taken from required surface mesh module attributes.

See also

PS 3.3 A.1.2.18 , A.57 and C.27

## 27.268.2 Member Enumeration Documentation

### 27.268.2.1 enum gdcm::Surface::STATES

Enumerator

***NO***

***YES***

***UNKNOWN***

***STATES\_END***

### 27.268.2.2 enum gdcm::Surface::VIEWType

Enumeration for Recommended Presentation [Type](#).

See also

Tag(0x0066, 0x000D) and PS 3.3 C.27.1.1.3

Enumerator

***SURFACE***

***WIREFRAME***

***POINTS***

***VIEWType\_END***

## 27.268.3 Constructor & Destructor Documentation

27.268.3.1 gdcm::Surface::Surface ( )

27.268.3.2 virtual gdcm::Surface::~~Surface ( ) [virtual]

## 27.268.4 Member Function Documentation

27.268.4.1 SegmentHelper::BasicCodedEntry const& gdcm::Surface::GetAlgorithmFamily ( ) const

27.268.4.2 SegmentHelper::BasicCodedEntry& gdcm::Surface::GetAlgorithmFamily ( )

27.268.4.3 const char\* gdcm::Surface::GetAlgorithmName ( ) const

27.268.4.4 const char\* gdcm::Surface::GetAlgorithmVersion ( ) const

27.268.4.5 const float\* gdcm::Surface::GetAxisOfRotation ( ) const

Note

Pointer is null if undefined

27.268.4.6 `const float* gdcm::Surface::GetCenterOfRotation ( ) const`

Note

Pointer is null if undefined

27.268.4.7 **STATES** `gdcm::Surface::GetFiniteVolume ( ) const`

27.268.4.8 **STATES** `gdcm::Surface::GetManifold ( ) const`

27.268.4.9 `float gdcm::Surface::GetMaximumPointDistance ( ) const`

27.268.4.10 `float gdcm::Surface::GetMeanPointDistance ( ) const`

27.268.4.11 **MeshPrimitive** `const& gdcm::Surface::GetMeshPrimitive ( ) const`

27.268.4.12 **MeshPrimitive** `& gdcm::Surface::GetMeshPrimitive ( )`

27.268.4.13 `unsigned long gdcm::Surface::GetNumberOfSurfacePoints ( ) const`

27.268.4.14 `unsigned long gdcm::Surface::GetNumberOfVectors ( ) const`

27.268.4.15 `const DataElement& gdcm::Surface::GetPointCoordinatesData ( ) const`

27.268.4.16 **DataElement** `& gdcm::Surface::GetPointCoordinatesData ( )`

27.268.4.17 `const float* gdcm::Surface::GetPointPositionAccuracy ( ) const`

Note

Pointer is null if undefined

27.268.4.18 `const float* gdcm::Surface::GetPointsBoundingBoxCoordinates ( ) const`

Note

Pointer is null if undefined

27.268.4.19 **SegmentHelper::BasicCodedEntry** `const& gdcm::Surface::GetProcessingAlgorithm ( ) const`

27.268.4.20 **SegmentHelper::BasicCodedEntry** `& gdcm::Surface::GetProcessingAlgorithm ( )`

27.268.4.21 `const unsigned short* gdcm::Surface::GetRecommendedDisplayCIELabValue ( ) const`

27.268.4.22 `unsigned short gdcm::Surface::GetRecommendedDisplayCIELabValue ( const unsigned int idx ) const`

27.268.4.23 `unsigned short gdcm::Surface::GetRecommendedDisplayGrayscaleValue ( ) const`

27.268.4.24 `float gdcm::Surface::GetRecommendedPresentationOpacity ( ) const`

- 27.268.4.25 **VIEWType** gdcm::Surface::GetRecommendedPresentationType ( ) const
- 27.268.4.26 **static STATES** gdcm::Surface::GetSTATES ( const char \* *state* ) [static]
- 27.268.4.27 **static const char\*** gdcm::Surface::GetSTATESString ( **STATES** *state* ) [static]
- 27.268.4.28 **const char\*** gdcm::Surface::GetSurfaceComments ( ) const
- 27.268.4.29 **unsigned long** gdcm::Surface::GetSurfaceNumber ( ) const
- 27.268.4.30 **bool** gdcm::Surface::GetSurfaceProcessing ( ) const
- 27.268.4.31 **const char\*** gdcm::Surface::GetSurfaceProcessingDescription ( ) const
- 27.268.4.32 **float** gdcm::Surface::GetSurfaceProcessingRatio ( ) const
- 27.268.4.33 **const float\*** gdcm::Surface::GetVectorAccuracy ( ) const
- 27.268.4.34 **const DataElement&** gdcm::Surface::GetVectorCoordinateData ( ) const
- 27.268.4.35 **DataElement&** gdcm::Surface::GetVectorCoordinateData ( )
- 27.268.4.36 **unsigned short** gdcm::Surface::GetVectorDimensionality ( ) const
- 27.268.4.37 **static VIEWType** gdcm::Surface::GetVIEWType ( const char \* *type* ) [static]
- 27.268.4.38 **static const char\*** gdcm::Surface::GetVIEWTypeString ( **VIEWType** *type* ) [static]
- 27.268.4.39 **void** gdcm::Surface::SetAlgorithmFamily ( **SegmentHelper::BasicCodedEntry** const & *BSE* )
- 27.268.4.40 **void** gdcm::Surface::SetAlgorithmName ( const char \* *str* )
- 27.268.4.41 **void** gdcm::Surface::SetAlgorithmVersion ( const char \* *str* )
- 27.268.4.42 **void** gdcm::Surface::SetAxisOfRotation ( const float \* *axis* )
- 27.268.4.43 **void** gdcm::Surface::SetCenterOfRotation ( const float \* *center* )
- 27.268.4.44 **void** gdcm::Surface::SetFiniteVolume ( **STATES** *state* )
- 27.268.4.45 **void** gdcm::Surface::SetManifold ( **STATES** *state* )
- 27.268.4.46 **void** gdcm::Surface::SetMaximumPointDistance ( float *maximum* )
- 27.268.4.47 **void** gdcm::Surface::SetMeanPointDistance ( float *average* )
- 27.268.4.48 **void** gdcm::Surface::SetMeshPrimitive ( **MeshPrimitive** & *mp* )
- 27.268.4.49 **void** gdcm::Surface::SetNumberOfSurfacePoints ( const unsigned long *nb* )
- 27.268.4.50 **void** gdcm::Surface::SetNumberOfVectors ( const unsigned long *nb* )

- 27.268.4.51 `void gdcM::Surface::SetPointCoordinatesData ( DataElement const & de )`
- 27.268.4.52 `void gdcM::Surface::SetPointPositionAccuracy ( const float * accuracies )`
- 27.268.4.53 `void gdcM::Surface::SetPointsBoundingBoxCoordinates ( const float * coordinates )`
- 27.268.4.54 `void gdcM::Surface::SetProcessingAlgorithm ( SegmentHelper::BasicCodedEntry const & BSE )`
- 27.268.4.55 `void gdcM::Surface::SetRecommendedDisplayCIELabValue ( const unsigned short vl[3] )`
- 27.268.4.56 `void gdcM::Surface::SetRecommendedDisplayCIELabValue ( const unsigned short vl, const unsigned int idx = 0 )`
- 27.268.4.57 `void gdcM::Surface::SetRecommendedDisplayCIELabValue ( const std::vector< unsigned short > & vl )`
- 27.268.4.58 `void gdcM::Surface::SetRecommendedDisplayGrayscaleValue ( const unsigned short vl )`
- 27.268.4.59 `void gdcM::Surface::SetRecommendedPresentationOpacity ( const float opacity )`
- 27.268.4.60 `void gdcM::Surface::SetRecommendedPresentationType ( VIEWType type )`
- 27.268.4.61 `void gdcM::Surface::SetSurfaceComments ( const char * comment )`
- 27.268.4.62 `void gdcM::Surface::SetSurfaceNumber ( const unsigned long nb )`
- 27.268.4.63 `void gdcM::Surface::SetSurfaceProcessing ( bool b )`
- 27.268.4.64 `void gdcM::Surface::SetSurfaceProcessingDescription ( const char * description )`
- 27.268.4.65 `void gdcM::Surface::SetSurfaceProcessingRatio ( const float ratio )`
- 27.268.4.66 `void gdcM::Surface::SetVectorAccuracy ( const float * accuracy )`
- 27.268.4.67 `void gdcM::Surface::SetVectorCoordinateData ( DataElement const & de )`
- 27.268.4.68 `void gdcM::Surface::SetVectorDimensionality ( const unsigned short dim )`

The documentation for this class was generated from the following file:

- [gdcMSurface.h](#)

## 27.269 gdcM::SurfaceHelper Class Reference

[SurfaceHelper](#) Helper class for [Surface](#) object.

```
#include <gdcMSurfaceHelper.h>
```

### Public Types

- `typedef std::vector< unsigned short > ColorArray`

## Static Public Member Functions

- `template<typename T , typename U >`  
`static std::vector< T > RecommendedDisplayCIELabToRGB (const ColorArray &CIELab, const U range←`  
`Max=255)`  
*Convert a DICOM CIE-Lab (after reading) color into RGB.*
- `template<typename U >`  
`static std::vector< float > RecommendedDisplayCIELabToRGB (const ColorArray &CIELab, const U range←`  
`Max=255)`  
*Convert a DICOM CIE-Lab (after reading) color into RGB.*
- `template<typename T , typename U >`  
`static ColorArray RGBToRecommendedDisplayCIELab (const std::vector< T > &RGB, const U rangeMax=255)`  
*Convert a RGB color into DICOM CIE-Lab (ready to write).*
- `template<typename T , typename U >`  
`static unsigned short RGBToRecommendedDisplayGrayscale (const std::vector< T > &RGB, const U range←`  
`Max=255)`  
*Convert a RGB color into DICOM grayscale (ready to write).*

### 27.269.1 Detailed Description

[SurfaceHelper](#) Helper class for [Surface](#) object.

### 27.269.2 Member Typedef Documentation

27.269.2.1 `typedef std::vector< unsigned short > gdcm::SurfaceHelper::ColorArray`

### 27.269.3 Member Function Documentation

27.269.3.1 `template<typename T , typename U > std::vector< T > gdcm::SurfaceHelper::RecommendedDisplayCIELabToRGB (`  
`const ColorArray & CIELab, const U rangeMax = 255 ) [static]`

Convert a DICOM CIE-Lab (after reading) color into RGB.

See also

PS 3.3 C.10.7.1.1

#### Parameters

<i>CIELab</i>	DICOM CIE-Lab array.
<i>rangeMax</i>	Max value of the RGB range.

#### Template Parameters

<i>T</i>	Type of CIELab components.
<i>U</i>	Type of rangeMax value.

27.269.3.2 `template<typename U > std::vector< float > gdcm::SurfaceHelper::RecommendedDisplayCIELabToRGB ( const`  
`ColorArray & CIELab, const U rangeMax = 255 ) [static]`

Convert a DICOM CIE-Lab (after reading) color into RGB.

See also

PS 3.3 C.10.7.1.1

Parameters

<i>CIELab</i>	DICOM CIE-Lab array.
<i>rangeMax</i>	Max value of the RGB range.

Template Parameters

<i>U</i>	Type of rangeMax value.
----------	-------------------------

```
27.269.3.3  template<typename T , typename U > SurfaceHelper::ColorArray gdcmm::SurfaceHelper::RGBToRecommendedDisplayCIELab ( const std::vector< T > & RGB, const U rangeMax = 255 )
[static]
```

Convert a RGB color into DICOM CIE-Lab (ready to write).

See also

PS 3.3 C.10.7.1.1

Parameters

<i>RGB</i>	RGB array.
<i>rangeMax</i>	Max value of the RGB range.

Template Parameters

<i>T</i>	Type of RGB components.
<i>U</i>	Type of rangeMax value.

```
27.269.3.4  template<typename T , typename U > unsigned short gdcmm::SurfaceHelper::RGBToRecommendedDisplayGrayscale (
const std::vector< T > & RGB, const U rangeMax = 255 ) [static]
```

Convert a RGB color into DICOM grayscale (ready to write).

See also

PS 3.3 C.27.1 tag(0062,000C)

Parameters

<i>RGB</i>	RGB array.
<i>rangeMax</i>	Max value of the RGB range.

Template Parameters

<i>T</i>	Type of RGB components.
----------	-------------------------



<i>U</i>	Type of rangeMax value.
----------	-------------------------

The documentation for this class was generated from the following file:


- [gdcmSurfaceHelper.h](#)

## 27.270 gdcm::SurfaceReader Class Reference


This class defines a SURFACE IE reader. It reads surface mesh module attributes.

```
#include <gdcmSurfaceReader.h>
```

Inheritance diagram for gdcm::SurfaceReader:

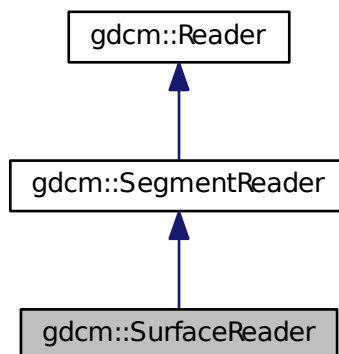

Collaboration diagram for `gdcM::SurfaceReader`:

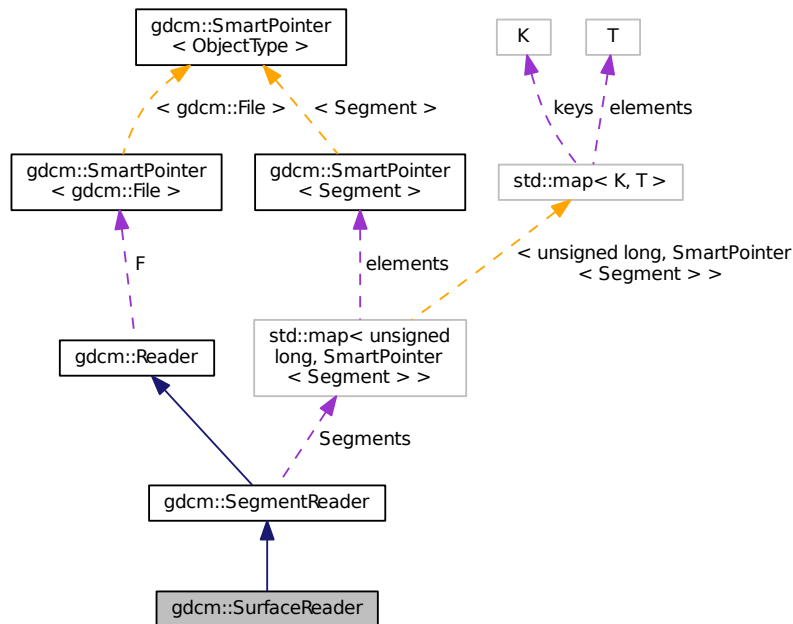

## Public Member Functions

- [SurfaceReader](#) ()
- virtual [~SurfaceReader](#) ()
- unsigned long [GetNumberOfSurfaces](#) () const
- virtual bool [Read](#) ()  
*Read.*

## Protected Member Functions

- bool [ReadPointMacro](#) ([SmartPointer< Surface >](#) surface, const [DataSet](#) &surfaceDS)
- bool [ReadSurface](#) (const [Item](#) &surfItem, const unsigned long idx)
- bool [ReadSurfaces](#) ()

## Additional Inherited Members

### 27.270.1 Detailed Description

This class defines a SURFACE IE reader. It reads surface mesh module attributes.

See also

PS 3.3 A.1.2.18 , A.57 and C.27

## 27.270.2 Constructor & Destructor Documentation

27.270.2.1 `gdcm::SurfaceReader::SurfaceReader ( )`

27.270.2.2 `virtual gdcm::SurfaceReader::~~SurfaceReader ( )` `[virtual]`

## 27.270.3 Member Function Documentation

27.270.3.1 `unsigned long gdcm::SurfaceReader::GetNumberOfSurfaces ( )` `const`

27.270.3.2 `virtual bool gdcm::SurfaceReader::Read ( )` `[virtual]`

Read.

Reimplemented from [gdcm::SegmentReader](#).

27.270.3.3 `bool gdcm::SurfaceReader::ReadPointMacro ( SmartPointer< Surface > surface, const DataSet & surfaceDS )`  
`[protected]`

27.270.3.4 `bool gdcm::SurfaceReader::ReadSurface ( const Item & surfaceItem, const unsigned long idx )` `[protected]`

27.270.3.5 `bool gdcm::SurfaceReader::ReadSurfaces ( )` `[protected]`

The documentation for this class was generated from the following file:

- [gdcmSurfaceReader.h](#)

## 27.271 gdcm::SurfaceWriter Class Reference

This class defines a SURFACE IE writer. It writes surface mesh module attributes.

```
#include <gdcmSurfaceWriter.h>
```

```
std::basic_ostream
```



- SurfaceWriter ()

- Write.*

- void ComputeNumberOf...

- unsigned long **NumberOfSurfaces**

## Additional Inherited Members

### 27.271.1 Detailed Description

This class defines a SURFACE IE writer. It writes surface mesh module attributes.

See also

PS 3.3 A.1.2.18 , A.57 and C.27

### 27.271.2 Constructor & Destructor Documentation

27.271.2.1 `gdcm::SurfaceWriter::SurfaceWriter ( )`

27.271.2.2 `virtual gdcm::SurfaceWriter::~~SurfaceWriter ( ) [virtual]`

### 27.271.3 Member Function Documentation

27.271.3.1 `void gdcm::SurfaceWriter::ComputeNumberOfSurfaces ( ) [protected]`

27.271.3.2 `unsigned long gdcm::SurfaceWriter::GetNumberOfSurfaces ( )`

27.271.3.3 `bool gdcm::SurfaceWriter::PrepareWrite ( ) [protected]`

27.271.3.4 `bool gdcm::SurfaceWriter::PrepareWritePointMacro ( SmartPointer< Surface > surface, DataSet & surfaceDS, const TransferSyntax & ts ) [protected]`

27.271.3.5 `void gdcm::SurfaceWriter::SetNumberOfSurfaces ( const unsigned long nb )`

27.271.3.6 `bool gdcm::SurfaceWriter::Write ( ) [virtual]`

Write.

Reimplemented from [gdcm::SegmentWriter](#).

### 27.271.4 Member Data Documentation

27.271.4.1 `unsigned long gdcm::SurfaceWriter::NumberOfSurfaces [protected]`

The documentation for this class was generated from the following file:

- [gdcmSurfaceWriter.h](#)

## 27.272 gdcm::SwapCode Class Reference

[SwapCode](#) representation.

```
#include <gdcmSwapCode.h>
```

## Public Types

- enum [SwapCodeType](#) {  
    [Unknown](#) = 0,  
    [LittleEndian](#) = 1234,  
    [BigEndian](#) = 4321,  
    [BadLittleEndian](#) = 3412,  
    [BadBigEndian](#) = 2143 }

## Public Member Functions

- [SwapCode](#) ([SwapCodeType](#) sc=[Unknown](#))
- [operator SwapCode::SwapCodeType](#) () const

## Static Public Member Functions

- static const char \* [GetSwapCodeString](#) ([SwapCode](#) const &sc)

## Static Protected Member Functions

- static int [GetIndex](#) ([SwapCode](#) const &sc)

## Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [SwapCode](#) &sc)

### 27.272.1 Detailed Description

[SwapCode](#) representation.

Examples:

[TestByteSwap.cxx](#).

### 27.272.2 Member Enumeration Documentation

#### 27.272.2.1 enum [gdcm::SwapCode::SwapCodeType](#)

Enumerator

***Unknown***

***LittleEndian***

***BigEndian***

***BadLittleEndian***

***BadBigEndian***

### 27.272.3 Constructor & Destructor Documentation

27.272.3.1 `gdcm::SwapCode::SwapCode ( SwapCodeType sc = Unknown )` `[inline]`

### 27.272.4 Member Function Documentation

27.272.4.1 `static int gdcm::SwapCode::GetIndex ( SwapCode const & sc )` `[static]`, `[protected]`

27.272.4.2 `static const char* gdcm::SwapCode::GetSwapCodeString ( SwapCode const & sc )` `[static]`

Referenced by `gdcm::operator<<()`.

27.272.4.3 `gdcm::SwapCode::operator SwapCode::SwapCodeType ( ) const` `[inline]`

### 27.272.5 Friends And Related Function Documentation

27.272.5.1 `std::ostream& operator<< ( std::ostream & os, const SwapCode & sc )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmSwapCode.h](#)

## 27.273 gdcm::SwapperDoOp Class Reference

```
#include <gdcmSwapper.h>
```

### Static Public Member Functions

- `template<typename T >`  
`static T Swap ( T val )`
- `template<typename T >`  
`static void SwapArray ( T *array, size_t n )`

### 27.273.1 Member Function Documentation

27.273.1.1 `template<typename T > static T gdcm::SwapperDoOp::Swap ( T val )` `[static]`

Referenced by `gdcm::Item::Read()`.

27.273.1.2 `template<typename T > static void gdcm::SwapperDoOp::SwapArray ( T * array, size_t n )` `[inline]`,  
`[static]`

The documentation for this class was generated from the following file:

- [gdcmSwapper.h](#)

## 27.274 gdcm::SwapperNoOp Class Reference

```
#include <gdcmSwapper.h>
```

### Static Public Member Functions

- template<typename T >  
static T [Swap](#) (T val)
- template<typename T >  
static void [SwapArray](#) (T \*, size\_t)

### 27.274.1 Detailed Description

Examples:

[ReadExplicitLengthSQIVR.cxx](#).

### 27.274.2 Member Function Documentation

**27.274.2.1** template<typename T > static T gdcm::SwapperNoOp::Swap ( T val ) [inline], [static]

Referenced by gdcm::EncodingImplementation< VR::VRBINARY >::Write().

**27.274.2.2** template<typename T > static void gdcm::SwapperNoOp::SwapArray ( T \*, size\_t ) [inline], [static]

Referenced by gdcm::EncodingImplementation< VR::VRBINARY >::Read().

The documentation for this class was generated from the following file:

- [gdcmSwapper.h](#)

## 27.275 gdcm::System Class Reference

Class to do system operation.

```
#include <gdcmSystem.h>
```

### Static Public Member Functions

- static bool [DeleteDirectory](#) (const char \*source)  
*remove a directory named source*
- static size\_t [EncodeBytes](#) (char \*out, const unsigned char \*data, int size)
- static bool [FileExists](#) (const char \*filename)  
*Check whether the specified file exist on the sytem.*
- static bool [FileIsDirectory](#) (const char \*name)  
*Check whether the file specified is a directory:*
- static bool [FileIsSymlink](#) (const char \*name)  
*Check whether name is a symlink.*



- static size\_t [FileSize](#) (const char \*filename)
- static time\_t [FileTime](#) (const char \*filename)
- static bool [FormatDateTime](#) (char date[22], time\_t t, long milliseconds=0)
- static bool [GetCurrentDateTime](#) (char date[22])
- static const char \* [GetCurrentModuleFileName](#) ()
- static const char \* [GetCurrentProcessFileName](#) ()
- static const char \* [GetCurrentResourcesDirectory](#) ()
- static const char \* [GetCWD](#) ()
- static bool [GetHostName](#) (char hostname[255])
- static const char \* [GetLastError](#) ()  
*Return the last error.*
- static const char \* [GetLocaleCharset](#) ()  
*return locale charmap*
- static const char \* [GetTimezoneOffsetFromUTC](#) ()
- static bool [MakeDirectory](#) (const char \*path)  
*Create a directory name path.*
- static bool [ParseDateTime](#) (time\_t &timep, const char date[22])  
*Parse a date stored as ASCII text into a time\_t structured (discard millisecond if any)*
- static bool [ParseDateTime](#) (time\_t &timep, long &milliseconds, const char date[22])
- static bool [RemoveFile](#) (const char \*source)  
*remove a file named source*
- static int [StrCaseCmp](#) (const char \*s1, const char \*s2)  
*consistent func for C99 spec of strcasecmp/strncasecmp*
- static int [StrNCaseCmp](#) (const char \*s1, const char \*s2, size\_t n)
- static char \* [StrSep](#) (char \*\*stringp, const char \*delim)
- static char \* [StrTokR](#) (char \*ptr, const char \*sep, char \*\*end)  
*strtok\_r*

## Static Protected Member Functions

- static bool [GetPermissions](#) (const char \*file, unsigned short &mode)  
*NOT THREAD SAFE.*
- static bool [SetPermissions](#) (const char \*file, unsigned short mode)

### 27.275.1 Detailed Description

Class to do system operation.

OS independent functionalities

### 27.275.2 Member Function Documentation

**27.275.2.1** static bool [gdcm::System::DeleteDirectory](#) ( const char \* *source* ) [static]

remove a directory named source

**27.275.2.2** static size\_t [gdcm::System::EncodeBytes](#) ( char \* *out*, const unsigned char \* *data*, int *size* ) [static]

Used internally by the [UIDGenerator](#) class to convert a uuid tape to a DICOM [VR:UI](#) type

**27.275.2.3** static bool gdcm::System::FileExists ( const char \* *filename* ) [static]

Check whether the specified file exist on the sytem.

Examples:

[EncapsulateFileInRawData.cxx](#), [gdcmorthoplanes.cxx](#), and [MagnifyFile.cxx](#).

**27.275.2.4** static bool gdcm::System::FileIsDirectory ( const char \* *name* ) [static]

Check whether the file specified is a directory:

Examples:

[gdcmorthoplanes.cxx](#), and [threadgdcm.cxx](#).

**27.275.2.5** static bool gdcm::System::FileIsSymlink ( const char \* *name* ) [static]

Check whether name is a symlink.

**27.275.2.6** static size\_t gdcm::System::FileSize ( const char \* *filename* ) [static]

Return the filesize. 0 if file does not exist.

Warning

you need to use FileExists to differentiate between empty file and missing file.  
for very large size file and on system where size\_t is not appropriate to store off\_t value the function will return 0.

Examples:

[CheckBigEndianBug.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [EncapsulateFileInRawData.cxx](#), and [SimpleScanner.cxx](#).

**27.275.2.7** static time\_t gdcm::System::FileTime ( const char \* *filename* ) [static]

Return the time of last modification of file 0 if the file does not exist

**27.275.2.8** static bool gdcm::System::FormatDateTime ( char *date[22]*, time\_t *t*, long *milliseconds* = 0 ) [static]

format as ASCII text a time\_t with milliseconds See [VR::DT](#) from DICOM PS 3.5 milliseconds is in the range [0, 999999]

**27.275.2.9** static bool gdcm::System::GetCurrentDateTime ( char *date[22]* ) [static]

Return the current data time, and format it as ASCII text. This is simply a call to gettimeofday + FormatDateTime, since WIN32 do not have an implementation for gettimeofday, this is more portable. The call time(0) is not precise for our resolution

27.275.2.10 `static const char* gdcm::System::GetCurrentModuleFileName ( ) [static]`

Return the directory the current module is located: NOT THREAD SAFE

27.275.2.11 `static const char* gdcm::System::GetCurrentProcessFileName ( ) [static]`

Return the directory the current process (executable) is located: NOT THREAD SAFE

27.275.2.12 `static const char* gdcm::System::GetCurrentResourcesDirectory ( ) [static]`

On some system (Apple) return the path to the current bundled 'Resources' directory NOT THREAD SAFE

27.275.2.13 `static const char* gdcm::System::GetCurrentWorkingDirectory ( ) [static]`

Return current working directory Warning: if current working path is too long (>2048 bytes) the call will fail and call will return NULL NOT THREAD SAFE

27.275.2.14 `static bool gdcm::System::GetHostName ( char hostname[255] ) [static]`

Retrieve the hostname, only the first 255 byte are copied. This may come handy to specify the Station Name

27.275.2.15 `static const char* gdcm::System::GetLastError ( ) [static]`

Return the last error.

27.275.2.16 `static const char* gdcm::System::GetLocaleCharSet ( ) [static]`

return locale charmap

27.275.2.17 `static bool gdcm::System::GetPermissions ( const char * file, unsigned short & mode ) [static],  
[protected]`

NOT THREAD SAFE.

27.275.2.18 `static const char* gdcm::System::GetTimezoneOffsetFromUTC ( ) [static]`

Return the value for Timezone Offset From UTC as string.

Warning

not thread safe

27.275.2.19 `static bool gdcm::System::MakeDirectory ( const char * path ) [static]`

Create a directory name path.

27.275.2.20 `static bool gdcmm::System::ParseDateTime ( time_t & timep, const char date[22] ) [static]`

Parse a date stored as ASCII text into a time\_t structured (discard millisecond if any)

27.275.2.21 `static bool gdcmm::System::ParseDateTime ( time_t & timep, long & milliseconds, const char date[22] ) [static]`

Parse a date stored as ASCII text into a time\_t structured and millisecond

See also

[FormatDateTime](#)

27.275.2.22 `static bool gdcmm::System::RemoveFile ( const char * source ) [static]`

remove a file named source

27.275.2.23 `static bool gdcmm::System::SetPermissions ( const char * file, unsigned short mode ) [static],  
[protected]`

27.275.2.24 `static int gdcmm::System::StrCaseCmp ( const char * s1, const char * s2 ) [static]`

consistent func for C99 spec of strcasecmp/strncasecmp

27.275.2.25 `static int gdcmm::System::StrNCaseCmp ( const char * s1, const char * s2, size_t n ) [static]`

Precondition

`n != 0`

27.275.2.26 `static char* gdcmm::System::StrSep ( char ** stringp, const char * delim ) [static]`

strsep param stringp is passed by pointer, it may be modified, you'll need to make a copy, in case you want to free the memory pointed at

27.275.2.27 `static char* gdcmm::System::StrTokR ( char * ptr, const char * sep, char ** end ) [static]`

strtok\_r

The documentation for this class was generated from the following file:

- [gdcmmSystem.h](#)

## 27.276 gdcmm::Table Class Reference

[Table.](#)

```
#include <gdcmmTable.h>
```

## Public Types

- typedef std::map< [Tag](#), [TableEntry](#) > [MapTableEntry](#)

## Public Member Functions

- [Table](#) ()
- [~Table](#) ()
- const [TableEntry](#) & [GetTableEntry](#) (const [Tag](#) &tag) const
- void [InsertEntry](#) ([Tag](#) const &tag, [TableEntry](#) const &te)

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Table](#) &\_val)

### 27.276.1 Detailed Description

[Table](#).

### 27.276.2 Member Typedef Documentation

27.276.2.1 typedef std::map<[Tag](#), [TableEntry](#)> [gdcm::Table::MapTableEntry](#)

### 27.276.3 Constructor & Destructor Documentation

27.276.3.1 [gdcm::Table::Table](#) ( ) [\[inline\]](#)

27.276.3.2 [gdcm::Table::~~Table](#) ( ) [\[inline\]](#)

### 27.276.4 Member Function Documentation

27.276.4.1 const [TableEntry](#)& [gdcm::Table::GetTableEntry](#) ( const [Tag](#) & *tag* ) const [\[inline\]](#)

27.276.4.2 void [gdcm::Table::InsertEntry](#) ( [Tag](#) const & *tag*, [TableEntry](#) const & *te* ) [\[inline\]](#)

### 27.276.5 Friends And Related Function Documentation

27.276.5.1 std::ostream& [operator<<](#) ( std::ostream &\_os, const [Table](#) &\_val ) [\[friend\]](#)

The documentation for this class was generated from the following file:

- [gdcmTable.h](#)

## 27.277 gdcm::TableEntry Class Reference

[TableEntry](#).

```
#include <gdcmTableEntry.h>
```

## Public Member Functions

- [TableEntry](#) (const char \*attribute=0, [Type](#) const &type=[Type](#)(), const char \*des=0)
- [~TableEntry](#) ()

### 27.277.1 Detailed Description

[TableEntry](#).

### 27.277.2 Constructor & Destructor Documentation

**27.277.2.1** `gdcm::TableEntry::TableEntry ( const char * attribute = 0, Type const & type = Type () , const char * des = 0 )`  
`[inline]`

**27.277.2.2** `gdcm::TableEntry::~~TableEntry ( )` `[inline]`

The documentation for this class was generated from the following file:


- [gdcmTableEntry.h](#)

## 27.278 gdcm::TableReader Class Reference

Class for representing a [TableReader](#).

```
#include <gdcmTableReader.h>
```

Inheritance diagram for `gdcm::TableReader`:

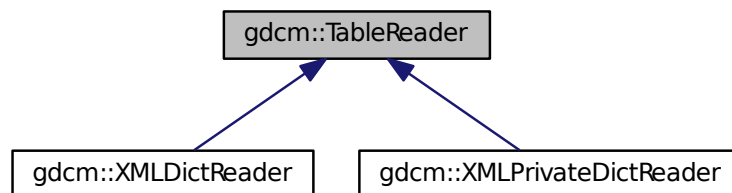

## Public Member Functions

- [TableReader](#) ([Defs](#) &defs)
- virtual [~TableReader](#) ()
- virtual void [CharacterDataHandler](#) (const char \*data, int length)
- virtual void [EndElement](#) (const char \*name)
- const [Defs](#) & [GetDefs](#) () const
- const char \* [GetFilename](#) ()

- void [HandleIOD](#) (const char \*\*atts)
- void [HandleIODEntry](#) (const char \*\*atts)
- void [HandleMacro](#) (const char \*\*atts)
- void [HandleMacroEntry](#) (const char \*\*atts)
- void [HandleMacroEntryDescription](#) (const char \*\*atts)
- void [HandleModule](#) (const char \*\*atts)
- void [HandleModuleEntry](#) (const char \*\*atts)
- void [HandleModuleEntryDescription](#) (const char \*\*atts)
- void [HandleModuleInclude](#) (const char \*\*atts)
- int [Read](#) ()
- void [SetFilename](#) (const char \*filename)
- virtual void [StartElement](#) (const char \*name, const char \*\*atts)

### 27.278.1 Detailed Description

Class for representing a [TableReader](#).

#### Note

This class is an empty shell meant to be derived

### 27.278.2 Constructor & Destructor Documentation

27.278.2.1 `gdcm::TableReader::TableReader ( Defs & defs )` `[inline]`

27.278.2.2 `virtual gdcm::TableReader::~~TableReader ( )` `[inline]`, `[virtual]`

### 27.278.3 Member Function Documentation

27.278.3.1 `virtual void gdcm::TableReader::CharacterDataHandler ( const char * data, int length )` `[virtual]`

Reimplemented in [gdcm::XMLDictReader](#), and [gdcm::XMLPrivateDictReader](#).

27.278.3.2 `virtual void gdcm::TableReader::EndElement ( const char * name )` `[virtual]`

Reimplemented in [gdcm::XMLDictReader](#), and [gdcm::XMLPrivateDictReader](#).

27.278.3.3 `const Defs& gdcm::TableReader::GetDefs ( ) const` `[inline]`

27.278.3.4 `const char* gdcm::TableReader::GetFilename ( )` `[inline]`

27.278.3.5 `void gdcm::TableReader::HandleIOD ( const char ** atts )`

27.278.3.6 `void gdcm::TableReader::HandleIODEntry ( const char ** atts )`

27.278.3.7 `void gdcm::TableReader::HandleMacro ( const char ** atts )`

27.278.3.8 `void gdcm::TableReader::HandleMacroEntry ( const char ** atts )`

```

27.278.3.9 void gdcM::TableReader::HandleMacroEntryDescription ( const char ** atts )
27.278.3.10 void gdcM::TableReader::HandleModule ( const char ** atts )
27.278.3.11 void gdcM::TableReader::HandleModuleEntry ( const char ** atts )
27.278.3.12 void gdcM::TableReader::HandleModuleEntryDescription ( const char ** atts )
27.278.3.13 void gdcM::TableReader::HandleModuleInclude ( const char ** atts )
27.278.3.14 int gdcM::TableReader::Read ( )
27.278.3.15 void gdcM::TableReader::SetFilename ( const char * filename ) [inline]
27.278.3.16 virtual void gdcM::TableReader::StartElement ( const char * name, const char ** atts ) [virtual]

```

Reimplemented in [gdcM::XMLDictReader](#), and [gdcM::XMLPrivateDictReader](#).


The documentation for this class was generated from the following file:

- [gdcMTableReader.h](#)

## 27.279 gdcM::network::TableRow Class Reference

```
#include <gdcMULTransitionTable.h>
```

Collaboration diagram for gdcM::network::TableRow:

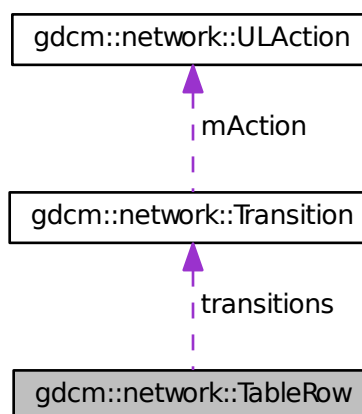

### Public Member Functions

- [TableRow](#) ()



- [~TableRow\(\)](#)

## Public Attributes

- [Transition](#) \* [transitions](#) [[cMaxStateID](#)]

## 27.279.1 Constructor & Destructor Documentation

27.279.1.1 `gdcM::network::TableRow::TableRow ( )` [[inline](#)]

References [gdcM::network::cMaxStateID](#), and [transitions](#).

27.279.1.2 `gdcM::network::TableRow::~~TableRow ( )` [[inline](#)]

References [gdcM::network::cMaxStateID](#), and [transitions](#).

## 27.279.2 Member Data Documentation

27.279.2.1 `Transition*` `gdcM::network::TableRow::transitions`[[cMaxStateID](#)]

Referenced by [TableRow\(\)](#), and [~TableRow\(\)](#).

The documentation for this class was generated from the following file:


- [gdcMULTransitionTable.h](#)

## 27.280 gdcM::Tag Class Reference

Class to represent a DICOM Data [Element](#) ([Attribute](#)) [Tag](#) (Group, [Element](#)). Basically an `uint32_t` which can also be expressed as two `uint16_t` (group and element)

```
#include <gdcMTag.h>
```

Inheritance diagram for `gdcM::Tag`:

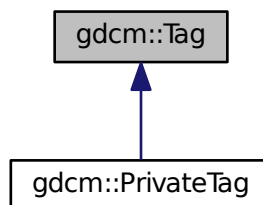

## Public Member Functions

- [Tag](#) (uint16\_t group, uint16\_t element)  
*Constructor with 2\*uint16\_t.*
- [Tag](#) (uint32\_t tag=0)  
*Constructor with 1\*uint32\_t Prefer the ctor that takes two uint16\_t.*
- [Tag](#) (const [Tag](#) &\_val)
- uint16\_t [GetElement](#) () const  
*Returns the 'Element number' of the given Tag.*
- uint32\_t [GetElementTag](#) () const  
*Returns the full tag value of the given Tag.*
- uint16\_t [GetGroup](#) () const  
*Returns the 'Group number' of the given Tag.*
- uint32\_t [GetLength](#) () const  
*return the length of tag (read: size on disk)*
- [Tag](#) [GetPrivateCreator](#) () const  
*Return the Private Creator Data Element tag of a private data element.*
- bool [IsGroupLength](#) () const  
*return whether the tag correspond to a group length tag:*
- bool [IsGroupXX](#) (const [Tag](#) &t) const  
*e.g 6002,3000 belong to groupXX: 6000,3000*
- bool [IsIllegal](#) () const  
*return if the tag is considered to be an illegal tag*
- bool [IsPrivate](#) () const
- bool [IsPrivateCreator](#) () const
- bool [IsPublic](#) () const
- bool [operator!=](#) (const [Tag](#) &\_val) const
- bool [operator<](#) (const [Tag](#) &\_val) const
- bool [operator<=](#) (const [Tag](#) &t2) const
- [Tag](#) & [operator=](#) (const [Tag](#) &\_val)
- bool [operator==](#) (const [Tag](#) &\_val) const
- const uint16\_t & [operator\[\]](#) (const unsigned int &\_id) const  
*Returns the Group or Element of the given Tag, depending on id (0/1)*
- uint16\_t & [operator\[\]](#) (const unsigned int &\_id)  
*Returns the Group or Element of the given Tag, depending on id (0/1)*
- std::string [PrintAsContinuousString](#) () const
- std::string [PrintAsContinuousUpperCaseString](#) () const  
*Same as PrintAsContinuousString, but hexadecimal [a-f] are printed using upper case.*
- std::string [PrintAsPipeSeparatedString](#) () const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)  
*Read a tag from binary representation.*
- bool [ReadFromCommaSeparatedString](#) (const char \*str)
- bool [ReadFromContinuousString](#) (const char \*str)
- bool [ReadFromPipeSeparatedString](#) (const char \*str)
- void [SetElement](#) (uint16\_t element)  
*Sets the 'Element number' of the given Tag.*
- void [SetElementTag](#) (uint16\_t group, uint16\_t element)

Sets the 'Group number' & 'Element number' of the given *Tag*.

- void [SetElementTag](#) (uint32\_t tag)

Sets the full tag value of the given *Tag*.

- void [SetGroup](#) (uint16\_t group)

Sets the 'Group number' of the given *Tag*.

- void [SetPrivateCreator](#) (*Tag* const &t)

Set private creator:

- template<typename TSwap >  
const std::ostream & [Write](#) (std::ostream &os) const

Write a tag in binary rep.

## Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const *Tag* &\_val)
- std::istream & [operator>>](#) (std::istream &\_is, *Tag* &\_val)

## 27.280.1 Detailed Description

Class to represent a DICOM Data [Element](#) ([Attribute](#)) *Tag* (Group, [Element](#)). Basically an uint32\_t which can also be expressed as two uint16\_t (group and element)

### Note

DATA ELEMENT TAG: A unique identifier for a Data [Element](#) composed of an ordered pair of numbers (a Group Number followed by an [Element](#) Number). GROUP NUMBER: The first number in the ordered pair of numbers that makes up a Data [Element Tag](#). ELEMENT NUMBER: The second number in the ordered pair of numbers that makes up a Data [Element Tag](#).

### Examples:

[ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJIPIDataSet.cxx](#), [DumpToSQLITE3.cxx](#), [DuplicatePCDE.cxx](#), [EncapsulateFileInRawData.cxx](#), [ExtractEncryptedContent.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [gdcmrtionplan.cxx](#), [gdcmrtplan.cxx](#), [GenAllIVR.cxx](#), [GenFakelIdentifyFile.cxx](#), [GenFakelImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetJPEGSamplePrecision.cxx](#), [GetSequenceUltrasound.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSExplicit.cxx](#), [MergeTwoFiles.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [PublicDict.cxx](#), [ReadAndDumpDICOMDIR.cxx](#), [ReadAndPrintAttributes.cxx](#), [ReadExplicitLengthSQIVR.cxx](#), [rle2img.cxx](#), [SimpleScanner.cxx](#), [SortImageReaderTest.cxx](#), [StreamImageReaderTest.cxx](#), [TraverseModules.cxx](#), and [VolumeSorter.cxx](#).

## 27.280.2 Constructor & Destructor Documentation

### 27.280.2.1 gdcm::Tag::Tag ( uint16\_t group, uint16\_t element ) [inline]

Constructor with 2\*uint16\_t.

### 27.280.2.2 gdcm::Tag::Tag ( uint32\_t tag = 0 ) [inline]

Constructor with 1\*uint32\_t Prefer the ctor that takes two uint16\_t.

**27.280.2.3** `gdcmm::Tag::Tag ( const Tag &_val ) [inline]`

References tag.

### 27.280.3 Member Function Documentation

**27.280.3.1** `uint16_t gdcmm::Tag::GetElement ( ) const [inline]`

Returns the 'Element number' of the given [Tag](#).

Examples:

[DuplicatePCDE.cxx](#), and [PublicDict.cxx](#).

Referenced by `gdcmm::DataSet::ComputeGroupLength()`, `IsGroupXX()`, `gdcmm::PrivateDict::PrintXML()`, `gdcmm::PrivateDict::PrivateTag()`, `gdcmm::SequenceOfFragments::ReadValue()`, and `SetPrivateCreator()`.

**27.280.3.2** `uint32_t gdcmm::Tag::GetElementTag ( ) const [inline]`

Returns the full tag value of the given [Tag](#).

**27.280.3.3** `uint16_t gdcmm::Tag::GetGroup ( ) const [inline]`

Returns the 'Group number' of the given [Tag](#).

Examples:

[DuplicatePCDE.cxx](#), and [GenAllVR.cxx](#).

Referenced by `gdcmm::DataSet::ComputeGroupLength()`, `gdcmm::CommandDataSet::Insert()`, `gdcmm::FileMetaInformation::Insert()`, `gdcmm::DataSet::Insert()`, `IsGroupXX()`, `gdcmm::PrivateDict::PrintXML()`, `gdcmm::SequenceOfFragments::ReadValue()`, `gdcmm::Attribute< Group, Element, TVR, TVM >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1 >::SetFromDataElement()`, `gdcmm::Attribute< Group, Element, TVR, VM::VM1_n >::SetFromDataElement()`, and `SetPrivateCreator()`.

**27.280.3.4** `uint32_t gdcmm::Tag::GetLength ( ) const [inline]`

return the length of tag (read: size on disk)

**27.280.3.5** `Tag gdcmm::Tag::GetPrivateCreator ( ) const [inline]`

Return the Private Creator Data [Element](#) tag of a private data element.

References `SetElement()`.

**27.280.3.6** `bool gdcmm::Tag::IsGroupLength ( ) const [inline]`

return whether the tag correspond to a group length tag:

**27.280.3.7** `bool gdcmm::Tag::IsGroupXX ( const Tag & t ) const` `[inline]`

e.g 6002,3000 belong to groupXX: 6000,3000

References `GetElement()`, `GetGroup()`, and `IsPrivate()`.

**27.280.3.8** `bool gdcmm::Tag::IsIllegal ( ) const` `[inline]`

return if the tag is considered to be an illegal tag

**27.280.3.9** `bool gdcmm::Tag::IsPrivate ( ) const` `[inline]`

PRIVATE DATA ELEMENT: Additional Data [Element](#), defined by an implementor, to communicate information that is not contained in Standard Data Elements. Private Data elements have odd Group Numbers.

Examples:

[DuplicatePCDE.cxx](#).

Referenced by `IsGroupXX()`, and `SetPrivateCreator()`.

**27.280.3.10** `bool gdcmm::Tag::IsPrivateCreator ( ) const` `[inline]`

Returns if tag is a Private Creator (xxxx,00yy), where xxxx is odd number and yy in [0x10,0xFF]

Examples:

[DuplicatePCDE.cxx](#).

**27.280.3.11** `bool gdcmm::Tag::IsPublic ( ) const` `[inline]`

STANDARD DATA ELEMENT: A Data [Element](#) defined in the DICOM Standard, and therefore listed in the DICOM Data [Element](#) Dictionary in PS 3.6. Is the [Tag](#) from the Public dict...well the implementation is buggy it does not prove the element is indeed in the dict...

**27.280.3.12** `bool gdcmm::Tag::operator!= ( const Tag & _val ) const` `[inline]`

References tag.

**27.280.3.13** `bool gdcmm::Tag::operator< ( const Tag & _val ) const` `[inline]`

DICOM Standard expects the Data [Element](#) to be sorted by Tags All other comparison can be constructed from this one and operator ==

References tag, and tags.

**27.280.3.14** `bool gdcmm::Tag::operator<= ( const Tag & t2 ) const` `[inline]`

**27.280.3.15** `Tag& gdcmm::Tag::operator= ( const Tag & _val )` `[inline]`

References tag.

27.280.3.16 `bool gdcmm::Tag::operator==( const Tag &_val ) const` `[inline]`

References tag.

27.280.3.17 `const uint16_t& gdcmm::Tag::operator[]( const unsigned int &_id ) const` `[inline]`

Returns the Group or [Element](#) of the given [Tag](#), depending on id (0/1)

27.280.3.18 `uint16_t& gdcmm::Tag::operator[]( const unsigned int &_id )` `[inline]`

Returns the Group or [Element](#) of the given [Tag](#), depending on id (0/1)

27.280.3.19 `std::string gdcmm::Tag::PrintAsContinuousString ( ) const`

Print tag value with no separating comma: eg. tag = "12345678" It comes in useful when reading tag values from XML file(in NativeDICOMModel)

27.280.3.20 `std::string gdcmm::Tag::PrintAsContinuousUpperCaseString ( ) const`

Same as `PrintAsContinuousString`, but hexadecimal [a-f] are printed using upper case.

27.280.3.21 `std::string gdcmm::Tag::PrintAsPipeSeparatedString ( ) const`

Print as a pipe separated string (GDCM 1.x compat only). Do not use in newer code

See also

[ReadFromPipeSeparatedString](#)

27.280.3.22 `template<typename TSwap> std::istream& gdcmm::Tag::Read ( std::istream &is )` `[inline]`

Read a tag from binary representation.

27.280.3.23 `bool gdcmm::Tag::ReadFromCommaSeparatedString ( const char * str )`

Read from a comma separated string. This is a highly user oriented function, the string should be formatted as↵ : 1234,5678 to specify the tag (0x1234,0x5678) The notation comes from the DICOM standard, and is handy to use from a command line program

27.280.3.24 `bool gdcmm::Tag::ReadFromContinuousString ( const char * str )`

Read From XML formatted tag value eg. tag = "12345678" It comes in useful when reading tag values from XML file(in NativeDICOMModel)

27.280.3.25 `bool gdcm::Tag::ReadFromPipeSeparatedString ( const char * str )`

Read from a pipe separated string (GDCM 1.x compat only). Do not use in newer code

See also

[ReadFromCommaSeparatedString](#)

27.280.3.26 `void gdcm::Tag::SetElement ( uint16_t element ) [inline]`

Sets the '[Element](#) number' of the given [Tag](#).

Examples:

[DuplicatePCDE.cxx](#), and [PublicDict.cxx](#).

Referenced by `GetPrivateCreator()`, and `gdcm::operator>>()`.

27.280.3.27 `void gdcm::Tag::SetElementTag ( uint16_t group, uint16_t element ) [inline]`

Sets the 'Group number' & '[Element](#) number' of the given [Tag](#).

27.280.3.28 `void gdcm::Tag::SetElementTag ( uint32_t tag ) [inline]`

Sets the full tag value of the given [Tag](#).

27.280.3.29 `void gdcm::Tag::SetGroup ( uint16_t group ) [inline]`

Sets the 'Group number' of the given [Tag](#).

Referenced by `gdcm::operator>>()`.

27.280.3.30 `void gdcm::Tag::SetPrivateCreator ( Tag const & t ) [inline]`

Set private creator:

Examples:

[DuplicatePCDE.cxx](#).

References `GetElement()`, `GetGroup()`, and `IsPrivate()`.

27.280.3.31 `template<typename TSwap> const std::ostream& gdcm::Tag::Write ( std::ostream & os ) const [inline]`

Write a tag in binary rep.

Referenced by `gdcm::SequenceOfItems::Write()`, `gdcm::Item::Write()`, and `gdcm::SequenceOfFragments::Write()`.

## 27.280.4 Friends And Related Function Documentation

27.280.4.1 `std::ostream& operator<< ( std::ostream &_os, const Tag &_val )` [[friend](#)]

27.280.4.2 `std::istream& operator>> ( std::istream &_is, Tag &_val )` [[friend](#)]

## 27.280.5 Member Data Documentation

27.280.5.1 `char gdcmm::Tag::bytes[4]`

27.280.5.2 `uint32_t gdcmm::Tag::tag`

Referenced by `operator!=()`, `operator<()`, `operator=()`, `operator==()`, and `Tag()`.

27.280.5.3 `uint16_t gdcmm::Tag::tags[2]`

Referenced by `operator<()`.

The documentation for this class was generated from the following file:

- [gdcmmTag.h](#)

## 27.281 gdcmm::TagPath Class Reference

class to handle a path of tag.

```
#include <gdcmmTagPath.h>
```

### Public Member Functions

- [TagPath](#) ()
- [~TagPath](#) ()
- bool [ConstructFromString](#) (const char \*path)
- bool [ConstructFromTagList](#) ([Tag](#) const \*l, unsigned int n)  
*Construct from a list of tags.*
- void [Print](#) (std::ostream &) const
- bool [Push](#) ([Tag](#) const &t)
- bool [Push](#) (unsigned int itemnum)

### Static Public Member Functions

- static bool [IsValid](#) (const char \*path)  
*Return if path is valid or not.*

### 27.281.1 Detailed Description

class to handle a path of tag.

Any Resemblance to Existing XPath is Purely Coincidental [ftp://medical.nema.org/medical/dicom/supps/sup118←\\_pc.pdf](ftp://medical.nema.org/medical/dicom/supps/sup118/_pc.pdf)



## 27.281.2 Constructor & Destructor Documentation

27.281.2.1 `gdcm::TagPath::TagPath ( )`

27.281.2.2 `gdcm::TagPath::~~TagPath ( )`

## 27.281.3 Member Function Documentation

27.281.3.1 `bool gdcm::TagPath::ConstructFromString ( const char * path )`

"/0018,0018/"... No space allowed, comma is use to separate tag group from tag element and slash is used to separate tag return false if invalid

27.281.3.2 `bool gdcm::TagPath::ConstructFromTagList ( Tag const * l, unsigned int n )`

Construct from a list of tags.

27.281.3.3 `static bool gdcm::TagPath::IsValid ( const char * path )` `[static]`

Return if path is valid or not.

27.281.3.4 `void gdcm::TagPath::Print ( std::ostream & ) const`

27.281.3.5 `bool gdcm::TagPath::Push ( Tag const & t )`

27.281.3.6 `bool gdcm::TagPath::Push ( unsigned int itemnum )`

The documentation for this class was generated from the following file:

- [gdcmTagPath.h](#)

## 27.282 gdcm::Testing Class Reference

class for testing

```
#include <gdcmTesting.h>
```

### Public Types

- `typedef const char *const (* MD5DataImagesType)[2]`
- `typedef const char *const (* MediaStorageDataFilesType)[2]`  
*return the table that map the media storage (as string) of a filename (gdcmData)*

### Public Member Functions

- `Testing ()`
- `~Testing ()`
- `void Print (std::ostream &os=std::cout)`  
*Print.*

## Static Public Member Functions

- static bool [ComputeFileMD5](#) (const char \*filename, char digest\_str[33])
- static bool [ComputeMD5](#) (const char \*buffer, unsigned long buf\_len, char digest\_str[33])
- static const char \* [GetDataExtraRoot](#) ()  
*Return the GDCM DATA EXTRA ROOT.*
- static const char \* [GetDataRoot](#) ()  
*Return the GDCM DATA ROOT.*
- static const char \* [GetFileName](#) (unsigned int file)
- static const char \*const \* [GetFileNames](#) ()  
*return the table of fullpath to gdcmData DICOM files:*
- static int [GetLossyFlagFromFile](#) (const char \*filepath)
- static const char \*const \* [GetMD5DataImage](#) (unsigned int file)
- static [MD5DataImagesType](#) [GetMD5DataImages](#) ()
- static const char \* [GetMD5FromBrokenFile](#) (const char \*filepath)
- static const char \* [GetMD5FromFile](#) (const char \*filepath)
- static const char \*const \* [GetMediaStorageDataFile](#) (unsigned int file)
- static [MediaStorageDataFilesType](#) [GetMediaStorageDataFiles](#) ()
- static const char \* [GetMediaStorageFromFile](#) (const char \*filepath)
- static unsigned int [GetNumberOfFileNames](#) ()
- static unsigned int [GetNumberOfMD5DataImages](#) ()
- static unsigned int [GetNumberOfMediaStorageDataFiles](#) ()
- static const char \* [GetPixelSpacingDataRoot](#) ()  
*Return the GDCM PIXEL SPACING DATA ROOT (See David Clunie website for dataset)*
- static std::streamoff [GetSelectedPrivateGroupOffsetFromFile](#) (const char \*filepath)
- static std::streamoff [GetSelectedTagsOffsetFromFile](#) (const char \*filepath)
- static const char \* [GetSourceDirectory](#) ()
- static std::streamoff [GetStreamOffsetFromFile](#) (const char \*filepath)
- static const char \* [GetTempDirectory](#) (const char \*subdir=0)
- static const wchar\_t \* [GetTempDirectoryW](#) (const wchar\_t \*subdir=0)  
*NOT THREAD SAFE.*
- static const char \* [GetTempFilename](#) (const char \*filename, const char \*subdir=0)  
*NOT THREAD SAFE.*
- static const wchar\_t \* [GetTempFilenameW](#) (const wchar\_t \*filename, const wchar\_t \*subdir=0)  
*NOT THREAD SAFE.*

### 27.282.1 Detailed Description

class for testing

this class is used for the nightly regression system for GDCM It makes heavily use of md5 computation

See also

[gdcm::MD5](#) class for md5 computation

### 27.282.2 Member Typedef Documentation

#### 27.282.2.1 typedef const char\* const(\* gdcm::Testing::MD5DataImagesType)[2]

return the table that map the md5 (as in md5sum) of the Pixel Data associated to a filename

27.282.2.2 `typedef const char* const(* gdcm::Testing::MediaStorageDataFileType)[2]`

return the table that map the media storage (as string) of a filename (gdcmData)

### 27.282.3 Constructor & Destructor Documentation

27.282.3.1 `gdcm::Testing::Testing ( ) [inline]`

27.282.3.2 `gdcm::Testing::~~Testing ( ) [inline]`

### 27.282.4 Member Function Documentation

27.282.4.1 `static bool gdcm::Testing::ComputeFileMD5 ( const char * filename, char digest_str[33] ) [static]`

Examples:

[MetalImageMD5Activiz.cs](#).

27.282.4.2 `static bool gdcm::Testing::ComputeMD5 ( const char * buffer, unsigned long buf_len, char digest_str[33] ) [static]`

MD5 stuff `digest_str` needs to be at least : `strlen = [2*16+1]`; string will be `\0` padded. (md5 are 32 bytes long) [Testing](#) is not meant to be shipped with an installed GDCM release, always prefer the [gdcm::MD5](#) API when doing md5 computation.

27.282.4.3 `static const char* gdcm::Testing::GetDataExtraRoot ( ) [static]`

Return the GDCM DATA EXTRA ROOT.

Examples:

[DiscriminateVolume.cxx](#), [reslicesphere.cxx](#), and [VolumeSorter.cxx](#).

27.282.4.4 `static const char* gdcm::Testing::GetDataRoot ( ) [static]`

Return the GDCM DATA ROOT.

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), and [MagnifyFile.cxx](#).

27.282.4.5 `static const char* gdcm::Testing::GetFileName ( unsigned int file ) [static]`

Examples:

[MetalImageMD5Activiz.cs](#).

27.282.4.6 static const char\* const\* gdcm::Testing::GetFileNames ( ) [static]

return the table of fullpath to gdcmData DICOM files:

Examples:

[TestReader.cxx](#).

27.282.4.7 static int gdcm::Testing::GetLossyFlagFromFile ( const char \* *filepath* ) [static]

Return the lossy flag of the given filename -1 -> Error 0 -> Lossless 1 -> Lossy

27.282.4.8 static const char\* const\* gdcm::Testing::GetMD5DataImage ( unsigned int *file* ) [static]

27.282.4.9 static MD5DataImagesType gdcm::Testing::GetMD5DataImages ( ) [static]

27.282.4.10 static const char\* gdcm::Testing::GetMD5FromBrokenFile ( const char \* *filepath* ) [static]

Return what should have been the md5 of file 'filepath' This is based on current GDCM implementation to decipher a broken DICOM file.

27.282.4.11 static const char\* gdcm::Testing::GetMD5FromFile ( const char \* *filepath* ) [static]

27.282.4.12 static const char\* const\* gdcm::Testing::GetMediaStorageDataFile ( unsigned int *file* ) [static]

27.282.4.13 static MediaStorageDataFilesType gdcm::Testing::GetMediaStorageDataFiles ( ) [static]

27.282.4.14 static const char\* gdcm::Testing::GetMediaStorageFromFile ( const char \* *filepath* ) [static]

Examples:

[MetaImageMD5Activiz.cs](#), and [TestReader.cxx](#).

27.282.4.15 static unsigned int gdcm::Testing::GetNumberOfFileNames ( ) [static]

Examples:

[MetaImageMD5Activiz.cs](#).

27.282.4.16 static unsigned int gdcm::Testing::GetNumberOfMD5DataImages ( ) [static]

27.282.4.17 static unsigned int gdcm::Testing::GetNumberOfMediaStorageDataFiles ( ) [static]

27.282.4.18 static const char\* gdcm::Testing::GetPixelSpacingDataRoot ( ) [static]

Return the GDCM PIXEL SPACING DATA ROOT (See David Clunie website for dataset)

27.282.4.19 `static std::streamoff gdcm::Testing::GetSelectedPrivateGroupOffsetFromFile ( const char * filepath ) [static]`

Return the offset just after private attribute (0009,0010,"GEMS\_IDEN\_01") if found. Otherwise the offset of the next attribute -1 if not found

27.282.4.20 `static std::streamoff gdcm::Testing::GetSelectedTagsOffsetFromFile ( const char * filepath ) [static]`

Return the offset just after Pixel Data Length (7fe0,0000) if found. Otherwise the offset of the very first pixel cell in Pixel Data -1 if not found

27.282.4.21 `static const char* gdcm::Testing::GetSourceDirectory ( ) [static]`

Examples:

[BasicAnonymizer.cs](#), and [ClinicalTrialIdentificationWorkflow.cs](#).

27.282.4.22 `static std::streamoff gdcm::Testing::GetStreamOffsetFromFile ( const char * filepath ) [static]`

Return the offset of the very first pixel cell in the PixelData -1 if not found

27.282.4.23 `static const char* gdcm::Testing::GetTempDirectory ( const char * subdir = 0 ) [static]`

NOT THREAD SAFE Returns the temp directory as used in testing needing to output data:

Examples:

[MetaImageMD5Activiz.cs](#).

27.282.4.24 `static const wchar_t* gdcm::Testing::GetTempDirectoryW ( const wchar_t * subdir = 0 ) [static]`

NOT THREAD SAFE.

27.282.4.25 `static const char* gdcm::Testing::GetTempFilename ( const char * filename, const char * subdir = 0 ) [static]`

NOT THREAD SAFE.

Examples:

[MetaImageMD5Activiz.cs](#).

27.282.4.26 `static const wchar_t* gdcm::Testing::GetTempFilenameW ( const wchar_t * filename, const wchar_t * subdir = 0 ) [static]`

NOT THREAD SAFE.

27.282.4.27 void `gdcM::Testing::Print ( std::ostream & os = std::cout )`

Print.

The documentation for this class was generated from the following file:

- [gdcMTesting.h](#)

## 27.283 gdcM::Trace Class Reference

[Trace](#).

```
#include <gdcMTrace.h>
```

### Public Member Functions

- [Trace](#) ()
- [~Trace](#) ()

### Static Public Member Functions

- static void [DebugOff](#) ()
- static void [DebugOn](#) ()
- static void [ErrorOff](#) ()
- static void [ErrorOn](#) ()
- static bool [GetDebugFlag](#) ()
- static std::ostream & [GetDebugStream](#) ()
- static bool [GetErrorFlag](#) ()
- static std::ostream & [GetErrorStream](#) ()
- static std::ostream & [GetStream](#) ()
- static bool [GetWarningFlag](#) ()
- static std::ostream & [GetWarningStream](#) ()
- static void [SetDebug](#) (bool debug)  
*Turn debug messages on (default: false)*
- static void [SetDebugStream](#) (std::ostream &os)  
*Explicitly set the stream which receive Debug messages:*
- static void [SetError](#) (bool debug)  
*Turn error messages on (default: true)*
- static void [SetErrorStream](#) (std::ostream &os)  
*Explicitly set the stream which receive Error messages:*
- static void [SetStream](#) (std::ostream &os)
- static void [SetStreamToFile](#) (const char \*filename)
- static void [SetWarning](#) (bool debug)  
*Turn warning messages on (default: true)*
- static void [SetWarningStream](#) (std::ostream &os)  
*Explicitly set the stream which receive Warning messages:*
- static void [WarningOff](#) ()
- static void [WarningOn](#) ()

### 27.283.1 Detailed Description

#### [Trace](#).

Debug / Warning and Error are encapsulated in this class by default the [Trace](#) class will redirect any debug/warning/error to `std::cerr`. Unless `SetStream` was specified with another (open) stream or `SetStreamToFile` was specified to a writable file on the system.

#### Warning

All string messages are removed during compilation time when compiled with `CMAKE_BUILD_TYPE` being set to either:

- Release
- MinSizeRel It is recommended to compile with `RelWithDebInfo` and/or `Debug` during prototyping of applications.

### 27.283.2 Constructor & Destructor Documentation

27.283.2.1 `gdcmm::Trace::Trace ( )`

27.283.2.2 `gdcmm::Trace::~~Trace ( )`

### 27.283.3 Member Function Documentation

27.283.3.1 `static void gdcmm::Trace::DebugOff ( )` [static]

Examples:

[MetalImageMD5Activiz.cs](#), and [TestReader.cxx](#).

27.283.3.2 `static void gdcmm::Trace::DebugOn ( )` [static]

Examples:

[Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

27.283.3.3 `static void gdcmm::Trace::ErrorOff ( )` [static]

Examples:

[MetalImageMD5Activiz.cs](#).

27.283.3.4 `static void gdcmm::Trace::ErrorOn ( )` [static]

27.283.3.5 `static bool gdcmm::Trace::GetDebugFlag ( )` [static]

27.283.3.6 `static std::ostream& gdcmm::Trace::GetDebugStream ( )` [static]

27.283.3.7 `static bool gdcmm::Trace::GetErrorFlag ( )` [static]

27.283.3.8 `static std::ostream& gdcm::Trace::GetErrorStream ( ) [static]`

27.283.3.9 `static std::ostream& gdcm::Trace::GetStream ( ) [static]`

27.283.3.10 `static bool gdcm::Trace::GetWarningFlag ( ) [static]`

27.283.3.11 `static std::ostream& gdcm::Trace::GetWarningStream ( ) [static]`

27.283.3.12 `static void gdcm::Trace::SetDebug ( bool debug ) [static]`

Turn debug messages on (default: false)

Examples:

[DumpToSQLITE3.cxx](#).

27.283.3.13 `static void gdcm::Trace::SetDebugStream ( std::ostream & os ) [static]`

Explicitely set the stream which receive Debug messages:

27.283.3.14 `static void gdcm::Trace::SetError ( bool debug ) [static]`

Turn error messages on (default: true)

27.283.3.15 `static void gdcm::Trace::SetErrorStream ( std::ostream & os ) [static]`

Explicitely set the stream which receive Error messages:

Examples:

[CStoreQtProgress.cxx](#).

27.283.3.16 `static void gdcm::Trace::SetStream ( std::ostream & os ) [static]`

Explicitely set the ostream for [gdcm::Trace](#) to report to This will set the DebugStream, WarningStream and ErrorStream at once:

27.283.3.17 `static void gdcm::Trace::SetStreamToFile ( const char * filename ) [static]`

Explicitely set the filename for [gdcm::Trace](#) to report to The file will be created (it will not append to existing file)

27.283.3.18 `static void gdcm::Trace::SetWarning ( bool debug ) [static]`

Turn warning messages on (default: true)

Examples:

[DumpToSQLITE3.cxx](#).



27.283.3.19 `static void gdcm::Trace::SetWarningStream ( std::ostream & os ) [static]`

Explicitely set the stream which receive Warning messages:

27.283.3.20 `static void gdcm::Trace::WarningOff ( ) [static]`

Examples:

[MetalImageMD5Activiz.cs](#), and [TestReader.cxx](#).

27.283.3.21 `static void gdcm::Trace::WarningOn ( ) [static]`

Examples:

[Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), and [StreamImageReaderTest.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmTrace.h](#)

## 27.284 gdcm::TransferSyntax Class Reference

Class to manipulate Transfer Syntax.

```
#include <gdcmTransferSyntax.h>
```

### Public Types

- enum [NegociatedType](#) {  
    [Unknown](#) = 0,  
    [Explicit](#),  
    [Implicit](#) }
- enum [TSType](#) {

```

ImplicitVRLittleEndian = 0,
ImplicitVRBigEndianPrivateGE,
ExplicitVRLittleEndian,
DeflatedExplicitVRLittleEndian,
ExplicitVRBigEndian,
JPEGBaselineProcess1,
JPEGExtendedProcess2_4,
JPEGExtendedProcess3_5,
JPEGsSpectralSelectionProcess6_8,
JPEGFullProgressionProcess10_12,
JPEGLosslessProcess14,
JPEGLosslessProcess14_1,
JPEGLSLossless,
JPEGLSNearLossless,
JPEG2000Lossless,
JPEG2000,
JPEG2000Part2Lossless,
JPEG2000Part2,
RLELossless,
MPEG2MainProfile,
ImplicitVRBigEndianACRNEMA,
CT_private_ELE,
JPIPReferenced,
TS_END }

```

## Public Member Functions

- [TransferSyntax](#) (TType type=[ImplicitVRLittleEndian](#))
- bool [CanStoreLossy](#) () const
- [NegociatedType](#) [GetNegociatedType](#) () const
- const char \* [GetString](#) () const
- [SwapCode](#) [GetSwapCode](#) () const
- bool [IsEncapsulated](#) () const
- bool [IsEncoded](#) () const
- bool [IsExplicit](#) () const
- bool [IsImplicit](#) () const
- bool [IsLossless](#) () const
- bool [IsLossy](#) () const
- bool [IsValid](#) () const
- [operator](#) TType () const

## Static Public Member Functions

- static const char \* [GetTSSString](#) (TType ts)
- static TType [GetTSType](#) (const char \*str)

## Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [TransferSyntax](#) &ts)

### 27.284.1 Detailed Description

Class to manipulate Transfer Syntax.

#### Note

TRANSFER SYNTAX (Standard and Private): A set of encoding rules that allow Application Entities to unambiguously negotiate the encoding techniques (e.g., Data [Element](#) structure, byte ordering, compression) they are able to support, thereby allowing these Application Entities to communicate.

**Todo** : The implementation is completely retarded -> see [gdcm::UIDs](#) for a replacement We need: IsSupported We need preprocess of raw/xml file We need GetFullName()

Need a notion of Private Syntax. As defined in PS 3.5. Section 9.2

See also

[UIDs](#)

Examples:

[GetJPEGSamplePrecision.cxx](#), and [LargeVRDSExplicit.cxx](#).

### 27.284.2 Member Enumeration Documentation

#### 27.284.2.1 enum gdcm::TransferSyntax::NegociatedType

Enumerator

***Unknown***

***Explicit***

***Implicit***

#### 27.284.2.2 enum gdcm::TransferSyntax::TSType

Enumerator

***ImplicitVRLittleEndian***

***ImplicitVRBigEndianPrivateGE***

***ExplicitVRLittleEndian***

***DeflatedExplicitVRLittleEndian***

***ExplicitVRBigEndian***

***JPEGBaselineProcess1***

***JPEGExtendedProcess2\_4***

***JPEGExtendedProcess3\_5***

***JPEGSpectralSelectionProcess6\_8***

***JPEGFullProgressionProcess10\_12***

***JPEGLosslessProcess14***

***JPEGLosslessProcess14\_1***

***JPEGLSLossless***

*JPEGLSNearLossless*  
*JPEG2000Lossless*  
*JPEG2000*  
*JPEG2000Part2Lossless*  
*JPEG2000Part2*  
*RLELossless*  
*MPEG2MainProfile*  
*ImplicitVRBigEndianACRNEMA*  
*CT\_private\_ELE*  
*JPIPReferenced*  
*TS\_END*

### 27.284.3 Constructor & Destructor Documentation

27.284.3.1 `gdcm::TransferSyntax::TransferSyntax ( TSType type = ImplicitVRLittleEndian )` `[inline]`

### 27.284.4 Member Function Documentation

27.284.4.1 `bool gdcm::TransferSyntax::CanStoreLossy ( ) const`

return true if TransFer Syntax Allow storing of Lossy Pixel Data

27.284.4.2 `NegotiatedType gdcm::TransferSyntax::GetNegociatedType ( ) const`

27.284.4.3 `const char* gdcm::TransferSyntax::GetString ( ) const` `[inline]`

References GetString().

27.284.4.4 `SwapCode gdcm::TransferSyntax::GetSwapCode ( ) const`

**Deprecated** Return the [SwapCode](#) associated with the Transfer Syntax. Be careful with the special GE private syntax the [DataSet](#) is written in little endian but the Pixel Data is in Big Endian.

27.284.4.5 `static const char* gdcm::TransferSyntax::GetTSString ( TSType ts )` `[static]`

Examples:

[LargeVRDSExplicit.cxx](#).

Referenced by GetString(), and gdcm::operator<<().

27.284.4.6 `static TSType gdcm::TransferSyntax::GetTSType ( const char * str )` `[static]`

27.284.4.7 `bool gdcm::TransferSyntax::IsEncapsulated ( ) const`

Examples:

[ExtractIconFromFile.cxx](#).

27.284.4.8 `bool gdcm::TransferSyntax::IsEncoded ( ) const`

27.284.4.9 `bool gdcm::TransferSyntax::IsExplicit ( ) const`

27.284.4.10 `bool gdcm::TransferSyntax::IsImplicit ( ) const`

27.284.4.11 `bool gdcm::TransferSyntax::IsLossless ( ) const`

Return true if the transfer syntax algorithm is a lossless algorithm

27.284.4.12 `bool gdcm::TransferSyntax::IsLossy ( ) const`

Return true if the transfer syntax algorithm is a lossy algorithm

27.284.4.13 `bool gdcm::TransferSyntax::IsValid ( ) const` `[inline]`

27.284.4.14 `gdcm::TransferSyntax::operator TSType ( ) const` `[inline]`

## 27.284.5 Friends And Related Function Documentation

27.284.5.1 `std::ostream& operator<< ( std::ostream & os, const TransferSyntax & ts )` `[friend]`

The documentation for this class was generated from the following file:

- [gdcmTransferSyntax.h](#)

## 27.285 gdcm::network::TransferSyntaxSub Class Reference

[TransferSyntaxSub](#) Table 9-15 TRANSFER SYNTAX SUB-ITEM FIELDS.

```
#include <gdcmTransferSyntaxSub.h>
```

### Public Member Functions

- [TransferSyntaxSub](#) ( )
- `const char * GetName ( ) const`
- `bool operator== (const TransferSyntaxSub &ts) const`
- `void Print (std::ostream &os) const`
- `std::istream & Read (std::istream &is)`
- `void SetName (const char *name)`
- `void SetNameFromUID (UIDs::TSName tsname)`
- `size_t Size ( ) const`
- `const std::ostream & Write (std::ostream &os) const`

### 27.285.1 Detailed Description

[TransferSyntaxSub](#) Table 9-15 TRANSFER SYNTAX SUB-ITEM FIELDS.

TODO what is the goal of :

[Table](#) 9-19 TRANSFER SYNTAX SUB-ITEM FIELDS

## 27.285.2 Constructor & Destructor Documentation

27.285.2.1 `gdcmm::network::TransferSyntaxSub::TransferSyntaxSub ( )`

## 27.285.3 Member Function Documentation

27.285.3.1 `const char* gdcmm::network::TransferSyntaxSub::GetName ( ) const` `[inline]`

27.285.3.2 `bool gdcmm::network::TransferSyntaxSub::operator== ( const TransferSyntaxSub & ts ) const` `[inline]`

27.285.3.3 `void gdcmm::network::TransferSyntaxSub::Print ( std::ostream & os ) const`

27.285.3.4 `std::istream& gdcmm::network::TransferSyntaxSub::Read ( std::istream & is )`

27.285.3.5 `void gdcmm::network::TransferSyntaxSub::SetName ( const char * name )`

27.285.3.6 `void gdcmm::network::TransferSyntaxSub::SetNameFromUID ( UIDs::TSName tsname )`

27.285.3.7 `size_t gdcmm::network::TransferSyntaxSub::Size ( ) const`

27.285.3.8 `const std::ostream& gdcmm::network::TransferSyntaxSub::Write ( std::ostream & os ) const`


The documentation for this class was generated from the following file:

- [gdcmmTransferSyntaxSub.h](#)

## 27.286 gdcmm::network::Transition Struct Reference

```
#include <gdcmmULTransitionTable.h>
```

Collaboration diagram for `gdcmm::network::Transition`:

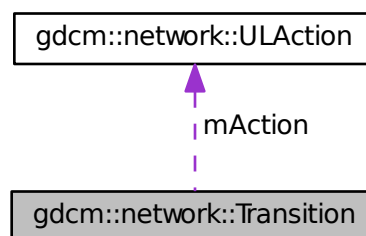

## Public Member Functions

- [Transition](#) ()

- [Transition](#) (int inEndState, [ULAction](#) \*inAction)
- [~Transition](#) ()

## Static Public Member Functions

- static [Transition](#) \* [MakeNew](#) (int inEndState, [ULAction](#) \*inAction)

## Public Attributes

- [ULAction](#) \* [mAction](#)
- int [mEnd](#)

## 27.286.1 Constructor & Destructor Documentation

27.286.1.1 `gdcm::network::Transition::Transition ( )` `[inline]`

References `gdcm::network::eStaDoesNotExist`.

Referenced by `MakeNew()`.

27.286.1.2 `gdcm::network::Transition::~~Transition ( )` `[inline]`

References `mAction`.

27.286.1.3 `gdcm::network::Transition::Transition ( int inEndState, ULAction * inAction )` `[inline]`

## 27.286.2 Member Function Documentation

27.286.2.1 `static Transition* gdcm::network::Transition::MakeNew ( int inEndState, ULAction * inAction )` `[inline]`,  
`[static]`

References `Transition()`.

## 27.286.3 Member Data Documentation

27.286.3.1 `ULAction* gdcm::network::Transition::mAction`

Referenced by `~Transition()`.

27.286.3.2 `int gdcm::network::Transition::mEnd`

The documentation for this struct was generated from the following file:

- [gdcmULTransitionTable.h](#)

## 27.287 `gdcm::Type` Class Reference

[Type.](#)

```
#include <gdcmType.h>
```

### Public Types

- enum [TypeType](#) {  
     [T1](#) = 0,  
     [T1C](#),  
     [T2](#),  
     [T2C](#),  
     [T3](#),  
     [UNKNOWN](#) }

### Public Member Functions

- [Type](#) ([TypeType](#) type=[UNKNOWN](#))
- [operator TypeType](#) () const

### Static Public Member Functions

- static const char \* [GetTypeString](#) ([TypeType](#) type)
- static [TypeType](#) [GetTypeType](#) (const char \*type)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [Type](#) &vr)

### 27.287.1 Detailed Description

[Type.](#)

#### Note

PS 3.5 7.4 DATA ELEMENT TYPE 7.4.1 TYPE 1 REQUIRED DATA ELEMENTS 7.4.2 TYPE 1C CONDITIONAL DATA ELEMENTS 7.4.3 TYPE 2 REQUIRED DATA ELEMENTS 7.4.4 TYPE 2C CONDITIONAL DATA ELEMENTS 7.4.5 TYPE 3 OPTIONAL DATA ELEMENTS

The intent of [Type](#) 2 Data Elements is to allow a zero length to be conveyed when the operator or application does not know its value or has a specific reason for not specifying its value. It is the intent that the device should support these Data Elements.

#### Examples:

[TraverseModules.cxx](#).



## 27.287.2 Member Enumeration Documentation

### 27.287.2.1 enum gdcm::Type::TypeType

Enumerator

***T1***  
***T1C***  
***T2***  
***T2C***  
***T3***  
***UNKNOWN***

## 27.287.3 Constructor & Destructor Documentation

27.287.3.1 gdcm::Type::Type ( TypeType *type* = UNKNOWN ) [inline]

## 27.287.4 Member Function Documentation

27.287.4.1 static const char\* gdcm::Type::GetTypeString ( TypeType *type* ) [static]

Referenced by gdcm::operator<<().

27.287.4.2 static TypeType gdcm::Type::GetTypeType ( const char \* *type* ) [static]

Referenced by gdcm::ModuleEntry::ModuleEntry().

27.287.4.3 gdcm::Type::operator TypeType ( ) const [inline]

## 27.287.5 Friends And Related Function Documentation

27.287.5.1 std::ostream& operator<< ( std::ostream & *os*, const Type & *vr* ) [friend]

The documentation for this class was generated from the following file:

- [gdcmType.h](#)

## 27.288 gdcm::UI Struct Reference

```
#include <gdcmVR.h>
```

### Public Attributes

- char [Internal](#) [64+1]

### Friends

- std::ostream & [operator<<](#) (std::ostream & *\_os*, const [UI](#) & *\_val*)

## 27.288.1 Friends And Related Function Documentation

27.288.1.1 `std::ostream& operator<< ( std::ostream &_os, const UI &_val )` [[friend](#)]

## 27.288.2 Member Data Documentation

27.288.2.1 `char gdcmm::UI::Internal[64+1]`

Referenced by `gdcmm::operator<<()`.

The documentation for this struct was generated from the following file:

- [gdcmmVR.h](#)

## 27.289 gdcmm::UIDGenerator Class Reference

Class for generating unique UID.

```
#include <gdcmmUIDGenerator.h>
```

### Public Member Functions

- [UIDGenerator](#) ()  
*By default the root of a UID is a GDCM Root...*
- `const char *` [Generate](#) ()

### Static Public Member Functions

- `static const char *` [GetGDCMUID](#) ()  
*Return the default (GDCM) root UID:*
- `static const char *` [GetRoot](#) ()
- `static bool` [IsValid](#) (`const char *uid`)
- `static void` [SetRoot](#) (`const char *root`)

### Static Protected Member Functions

- `static bool` [GenerateUUID](#) (`unsigned char *uuid_data`)

## 27.289.1 Detailed Description

Class for generating unique UID.

## Note

bla **Usage:** When constructing a [Series](#) or [Study](#) UID, user *has* to keep around the UID, otherwise the UID Generator will simply forget the value and create a new UID.

## Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenFakeImage.cxx](#), [GetSubSequenceData.cxx](#), [StreamImageReaderTest.cxx](#), and [uid\\_unique.cxx](#).

## 27.289.2 Constructor &amp; Destructor Documentation

27.289.2.1 `gdcm::UIDGenerator::UIDGenerator ( )` `[inline]`

By default the root of a UID is a GDCM Root...

## 27.289.3 Member Function Documentation

27.289.3.1 `const char* gdcm::UIDGenerator::Generate ( )`

Internally uses a `std::string`, so two calls have the same pointer ! save into a `std::string` In summary do not write code like that: `const char *uid1 = uid.Generate(); const char *uid2 = uid.Generate();` since `uid1 == uid2`

## Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GenFakeImage.cxx](#), [GetSubSequenceData.cxx](#), [StreamImageReaderTest.cxx](#), and [uid\\_unique.cxx](#).

27.289.3.2 `static bool gdcm::UIDGenerator::GenerateUUID ( unsigned char * uuid_data )` `[static]`, `[protected]`

27.289.3.3 `static const char* gdcm::UIDGenerator::GetGDCMUID ( )` `[static]`

Return the default (GDCM) root UID:

27.289.3.4 `static const char* gdcm::UIDGenerator::GetRoot ( )` `[static]`

## Examples:

[ClinicalTrialIdentificationWorkflow.cs](#), [ReformatFile.cs](#), and [StandardizeFiles.cs](#).

27.289.3.5 `static bool gdcm::UIDGenerator::IsValid ( const char * uid )` `[static]`

Find out if the string is a valid UID or not

**Todo** : Move that in `DataStructureAndEncoding` (see `FileMetaInformation::CheckFileMetaInformation`)

27.289.3.6 static void gdcm::UIDGenerator::SetRoot ( const char \* root ) [static]

The current implementation in GDCM make use of the UUID implementation (RFC 4122) and has been successfully been tested for a root of size 26 bytes. Any longer root should work (the [Generate\(\)](#) function will return a string), but will truncate the high bits of the 128bits UUID until the generated string fits on 64 bits. The authors disclaims any responsibility for guaranteeing uniqueness of [UIDs](#) when the root is longer than 26 bytes.

#### Examples:

[ClinicalTrialIdentificationWorkflow.cs](#), [ReformatFile.cs](#), [StandardizeFiles.cs](#), and [uid\\_unique.cxx](#).

The documentation for this class was generated from the following file:

- [gdcmUIDGenerator.h](#)

## 27.290 gdcm::UIDs Class Reference

all known uids

```
#include <gdcmUIDs.h>
```

### Public Types

- typedef const char \*const (\* [TransferSyntaxStringsType](#))[2]
- enum [TSName](#) {  
[VerificationSOPClass](#) = 1,  
[ImplicitVRLittleEndianDefaultTransferSyntaxforDICOM](#) = 2,  
[ExplicitVRLittleEndian](#) = 3,  
[DeflatedExplicitVRLittleEndian](#) = 4,  
[ExplicitVRBigEndian](#) = 5,  
[JPEGBaselineProcess1DefaultTransferSyntaxforLossyJPEG8BitImageCompression](#) = 6,  
[JPEGExtendedProcess24DefaultTransferSyntaxforLossyJPEG12BitImageCompressionProcess4only](#) = 7,  
[JPEGExtendedProcess35Retired](#) = 8,  
[JPEGSpectralSelectionNonHierarchicalProcess68Retired](#) = 9,  
[JPEGSpectralSelectionNonHierarchicalProcess79Retired](#) = 10,  
[JPEGFullProgressionNonHierarchicalProcess1012Retired](#) = 11,  
[JPEGFullProgressionNonHierarchicalProcess1113Retired](#) = 12,  
[JPEGLosslessNonHierarchicalProcess14](#) = 13,  
[JPEGLosslessNonHierarchicalProcess15Retired](#) = 14,  
[JPEGExtendedHierarchicalProcess1618Retired](#) = 15,  
[JPEGExtendedHierarchicalProcess1719Retired](#) = 16,  
[JPEGSpectralSelectionHierarchicalProcess2022Retired](#) = 17,  
[JPEGSpectralSelectionHierarchicalProcess2123Retired](#) = 18,  
[JPEGFullProgressionHierarchicalProcess2426Retired](#) = 19,  
[JPEGFullProgressionHierarchicalProcess2527Retired](#) = 20,  
[JPEGLosslessHierarchicalProcess28Retired](#) = 21,  
[JPEGLosslessHierarchicalProcess29Retired](#) = 22,  
[JPEGLosslessNonHierarchicalFirstOrderPredictionProcess14SelectionValue1DefaultTransferSyntaxforLossless](#)↵

JPEGImageCompression = 23,  
JPEGLSLosslessImageCompression = 24,  
JPEGLSLossyNearLosslessImageCompression = 25,  
JPEG2000ImageCompressionLosslessOnly = 26,  
JPEG2000ImageCompression = 27,  
JPEG2000Part2MulticomponentImageCompressionLosslessOnly = 28,  
JPEG2000Part2MulticomponentImageCompression = 29,  
JPIPReferenced = 30,  
JPIPReferencedDeflate = 31,  
MPEG2MainProfileMainLevel = 32,  
RLELossless = 33,  
RFC2557MIMEencapsulation = 34,  
XMLEncoding = 35,  
MediaStorageDirectoryStorage = 36,  
TalairachBrainAtlasFrameofReference = 37,  
SPM2T1FrameofReference = 38,  
SPM2T2FrameofReference = 39,  
SPM2PDFFrameofReference = 40,  
SPM2EPIFrameofReference = 41,  
SPM2FILT1FrameofReference = 42,  
SPM2PETFrameofReference = 43,  
SPM2TRANSMFrameofReference = 44,  
SPM2SPECTFrameofReference = 45,  
SPM2GRAYFrameofReference = 46,  
SPM2WHITEFrameofReference = 47,  
SPM2CSFFFrameofReference = 48,  
SPM2BRAINMASKFrameofReference = 49,  
SPM2AVG305T1FrameofReference = 50,  
SPM2AVG152T1FrameofReference = 51,  
SPM2AVG152T2FrameofReference = 52,  
SPM2AVG152PDFFrameofReference = 53,  
SPM2SINGLESUBJT1FrameofReference = 54,  
ICBM452T1FrameofReference = 55,  
ICBMSingleSubjectMRIFrameofReference = 56,  
BasicStudyContentNotificationSOPClassRetired = 57,  
StorageCommitmentPushModelSOPClass = 58,  
StorageCommitmentPushModelSOPInstance = 59,  
StorageCommitmentPullModelSOPClassRetired = 60,  
StorageCommitmentPullModelSOPInstanceRetired = 61,  
ProceduralEventLoggingSOPClass = 62,  
ProceduralEventLoggingSOPInstance = 63,  
SubstanceAdministrationLoggingSOPClass = 64,  
SubstanceAdministrationLoggingSOPInstance = 65,  
DICOMUIDRegistry = 66,  
DICOMControlledTerminology = 67,  
DICOMApplicationContextName = 68,  
DetachedPatientManagementSOPClassRetired = 69,  
DetachedPatientManagementMetaSOPClassRetired = 70,  
DetachedVisitManagementSOPClassRetired = 71,  
DetachedStudyManagementSOPClassRetired = 72,  
StudyComponentManagementSOPClassRetired = 73,  
ModalityPerformedProcedureStepSOPClass = 74,  
ModalityPerformedProcedureStepRetrieveSOPClass = 75,  
ModalityPerformedProcedureStepNotificationSOPClass = 76,  
DetachedResultsManagementSOPClassRetired = 77,  
DetachedResultsManagementMetaSOPClassRetired = 78,  
DetachedStudyManagementMetaSOPClassRetired = 79,  
DetachedInterpretationManagementSOPClassRetired = 80,  
StorageServiceClass = 81,  
BasicFilmSessionSOPClass = 82,

[BreastTomosynthesisImageStorage](#) }

• enum [TSType](#) {

```
uid_1_2_840_10008_1_1 = 1,  
uid_1_2_840_10008_1_2 = 2,  
uid_1_2_840_10008_1_2_1 = 3,  
uid_1_2_840_10008_1_2_1_99 = 4,  
uid_1_2_840_10008_1_2_2 = 5,  
uid_1_2_840_10008_1_2_4_50 = 6,  
uid_1_2_840_10008_1_2_4_51 = 7,  
uid_1_2_840_10008_1_2_4_52 = 8,  
uid_1_2_840_10008_1_2_4_53 = 9,  
uid_1_2_840_10008_1_2_4_54 = 10,  
uid_1_2_840_10008_1_2_4_55 = 11,  
uid_1_2_840_10008_1_2_4_56 = 12,  
uid_1_2_840_10008_1_2_4_57 = 13,  
uid_1_2_840_10008_1_2_4_58 = 14,  
uid_1_2_840_10008_1_2_4_59 = 15,  
uid_1_2_840_10008_1_2_4_60 = 16,  
uid_1_2_840_10008_1_2_4_61 = 17,  
uid_1_2_840_10008_1_2_4_62 = 18,  
uid_1_2_840_10008_1_2_4_63 = 19,  
uid_1_2_840_10008_1_2_4_64 = 20,  
uid_1_2_840_10008_1_2_4_65 = 21,  
uid_1_2_840_10008_1_2_4_66 = 22,  
uid_1_2_840_10008_1_2_4_70 = 23,  
uid_1_2_840_10008_1_2_4_80 = 24,  
uid_1_2_840_10008_1_2_4_81 = 25,  
uid_1_2_840_10008_1_2_4_90 = 26,  
uid_1_2_840_10008_1_2_4_91 = 27,  
uid_1_2_840_10008_1_2_4_92 = 28,  
uid_1_2_840_10008_1_2_4_93 = 29,  
uid_1_2_840_10008_1_2_4_94 = 30,  
uid_1_2_840_10008_1_2_4_95 = 31,  
uid_1_2_840_10008_1_2_4_100 = 32,  
uid_1_2_840_10008_1_2_5 = 33,  
uid_1_2_840_10008_1_2_6_1 = 34,  
uid_1_2_840_10008_1_2_6_2 = 35,  
uid_1_2_840_10008_1_3_10 = 36,  
uid_1_2_840_10008_1_4_1_1 = 37,  
uid_1_2_840_10008_1_4_1_2 = 38,  
uid_1_2_840_10008_1_4_1_3 = 39,  
uid_1_2_840_10008_1_4_1_4 = 40,  
uid_1_2_840_10008_1_4_1_5 = 41,  
uid_1_2_840_10008_1_4_1_6 = 42,  
uid_1_2_840_10008_1_4_1_7 = 43,  
uid_1_2_840_10008_1_4_1_8 = 44,  
uid_1_2_840_10008_1_4_1_9 = 45,  
uid_1_2_840_10008_1_4_1_10 = 46,  
uid_1_2_840_10008_1_4_1_11 = 47,  
uid_1_2_840_10008_1_4_1_12 = 48,  
uid_1_2_840_10008_1_4_1_13 = 49,  
uid_1_2_840_10008_1_4_1_14 = 50,  
uid_1_2_840_10008_1_4_1_15 = 51,  
uid_1_2_840_10008_1_4_1_16 = 52,  
uid_1_2_840_10008_1_4_1_17 = 53,  
uid_1_2_840_10008_1_4_1_18 = 54,  
uid_1_2_840_10008_1_4_2_1 = 55,  
uid_1_2_840_10008_1_4_2_2 = 56,  
uid_1_2_840_10008_1_9 = 57,  
uid_1_2_840_10008_1_20_1 = 58,  
uid_1_2_840_10008_1_20_1_1 = 59,  
uid_1_2_840_10008_1_20_2 = 60,
```

```
uid_1_2_840_10008_5_1_4_1_1_13_1_3 }
```

## Public Member Functions

- const char \* [GetName](#) () const
- const char \* [GetString](#) () const
- [operator TSType](#) () const
- bool [SetFromUID](#) (const char \*str)

## Static Public Member Functions

- static unsigned int [GetNumberOfTransferSyntaxStrings](#) ()
- static const char \*const \* [GetTransferSyntaxString](#) (unsigned int ts)
- static [TransferSyntaxStringsType](#) [GetTransferSyntaxStrings](#) ()
- static const char \* [GetUIDName](#) (unsigned int ts)
- static const char \* [GetUIDString](#) (unsigned int ts)

### 27.290.1 Detailed Description

all known uids

Examples:

[GenerateStandardSOPClasses.cxx](#).

### 27.290.2 Member Typedef Documentation

27.290.2.1 `typedef const char* const(* gdcmm::UIDs::TransferSyntaxStringsType)[2]`

### 27.290.3 Member Enumeration Documentation

27.290.3.1 `enum gdcmm::UIDs::TSName`

Enumerator

***VerificationSOPClass***

***ImplicitVRLittleEndianDefaultTransferSyntaxforDICOM***

***ExplicitVRLittleEndian***

***DeflatedExplicitVRLittleEndian***

***ExplicitVRBigEndian***

***JPEGBaselineProcess1DefaultTransferSyntaxforLossyJPEG8BitImageCompression***

***JPEGExtendedProcess24DefaultTransferSyntaxforLossyJPEG12BitImageCompressionProcess4only***

***JPEGExtendedProcess35Retired***

***JPEGSpectralSelectionNonHierarchicalProcess68Retired***

***JPEGSpectralSelectionNonHierarchicalProcess79Retired***

***JPEGFullProgressionNonHierarchicalProcess1012Retired***

***JPEGFullProgressionNonHierarchicalProcess1113Retired***



*JPEGLosslessNonHierarchicalProcess14*  
*JPEGLosslessNonHierarchicalProcess15Retired*  
*JPEGExtendedHierarchicalProcess1618Retired*  
*JPEGExtendedHierarchicalProcess1719Retired*  
*JPEGSpectralSelectionHierarchicalProcess2022Retired*  
*JPEGSpectralSelectionHierarchicalProcess2123Retired*  
*JPEGFullProgressionHierarchicalProcess2426Retired*  
*JPEGFullProgressionHierarchicalProcess2527Retired*  
*JPEGLosslessHierarchicalProcess28Retired*  
*JPEGLosslessHierarchicalProcess29Retired*  
*JPEGLosslessNonHierarchicalFirstOrderPredictionProcess14SelectionValue1DefaultTransferSyntaxforLosslessJPEGImage*

*JPEGLSLosslessImageCompression*  
*JPEGLSLossyNearLosslessImageCompression*  
*JPEG2000ImageCompressionLosslessOnly*  
*JPEG2000ImageCompression*  
*JPEG2000Part2MulticomponentImageCompressionLosslessOnly*  
*JPEG2000Part2MulticomponentImageCompression*  
*JPIPRreferenced*  
*JPIPRreferencedDeflate*  
*MPEG2MainProfileMainLevel*  
*RLELossless*  
*RFC2557MIMEencapsulation*  
*XMLEncoding*  
*MediaStorageDirectoryStorage*  
*TalairachBrainAtlasFrameofReference*  
*SPM2T1FrameofReference*  
*SPM2T2FrameofReference*  
*SPM2PDFFrameofReference*  
*SPM2EPIFrameofReference*  
*SPM2FIL T1FrameofReference*  
*SPM2PETFrameofReference*  
*SPM2TRANSMFrameofReference*  
*SPM2SPECTFrameofReference*  
*SPM2GRAYFrameofReference*  
*SPM2WHITEFrameofReference*  
*SPM2CSFFFrameofReference*  
*SPM2BRAINMASKFrameofReference*  
*SPM2AVG305T1FrameofReference*  
*SPM2AVG152T1FrameofReference*  
*SPM2AVG152T2FrameofReference*  
*SPM2AVG152PDFFrameofReference*

*SPM2SINGLESUBJT1FrameofReference*  
*ICBM452T1FrameofReference*  
*ICBMSingleSubjectMRIFrameofReference*  
*BasicStudyContentNotificationSOPClassRetired*  
*StorageCommitmentPushModelSOPClass*  
*StorageCommitmentPushModelSOPInstance*  
*StorageCommitmentPullModelSOPClassRetired*  
*StorageCommitmentPullModelSOPInstanceRetired*  
*ProceduralEventLoggingSOPClass*  
*ProceduralEventLoggingSOPInstance*  
*SubstanceAdministrationLoggingSOPClass*  
*SubstanceAdministrationLoggingSOPInstance*  
*DICOMUIDRegistry*  
*DICOMControlledTerminology*  
*DICOMApplicationContextName*  
*DetachedPatientManagementSOPClassRetired*  
*DetachedPatientManagementMetaSOPClassRetired*  
*DetachedVisitManagementSOPClassRetired*  
*DetachedStudyManagementSOPClassRetired*  
*StudyComponentManagementSOPClassRetired*  
*ModalityPerformedProcedureStepSOPClass*  
*ModalityPerformedProcedureStepRetrieveSOPClass*  
*ModalityPerformedProcedureStepNotificationSOPClass*  
*DetachedResultsManagementSOPClassRetired*  
*DetachedResultsManagementMetaSOPClassRetired*  
*DetachedStudyManagementMetaSOPClassRetired*  
*DetachedInterpretationManagementSOPClassRetired*  
*StorageServiceClass*  
*BasicFilmSessionSOPClass*  
*BasicFilmBoxSOPClass*  
*BasicGrayscaleImageBoxSOPClass*  
*BasicColorImageBoxSOPClass*  
*ReferencedImageBoxSOPClassRetired*  
*BasicGrayscalePrintManagementMetaSOPClass*  
*ReferencedGrayscalePrintManagementMetaSOPClassRetired*  
*PrintJobSOPClass*  
*BasicAnnotationBoxSOPClass*  
*PrinterSOPClass*  
*PrinterConfigurationRetrievalSOPClass*  
*PrinterSOPInstance*  
*PrinterConfigurationRetrievalSOPInstance*  
*BasicColorPrintManagementMetaSOPClass*

*ReferencedColorPrintManagementMetaSOPClassRetired*  
*VOILUTBoxSOPClass*  
*PresentationLUTSOPClass*  
*ImageOverlayBoxSOPClassRetired*  
*BasicPrintImageOverlayBoxSOPClassRetired*  
*PrintQueueSOPInstanceRetired*  
*PrintQueueManagementSOPClassRetired*  
*StoredPrintStorageSOPClassRetired*  
*HardcopyGrayscaleImageStorageSOPClassRetired*  
*HardcopyColorImageStorageSOPClassRetired*  
*PullPrintRequestSOPClassRetired*  
*PullStoredPrintManagementMetaSOPClassRetired*  
*MediaCreationManagementSOPClassUID*  
*ComputedRadiographyImageStorage*  
*DigitalXRayImageStorageForPresentation*  
*DigitalXRayImageStorageForProcessing*  
*DigitalMammographyXRayImageStorageForPresentation*  
*DigitalMammographyXRayImageStorageForProcessing*  
*DigitalIntraoralXRayImageStorageForPresentation*  
*DigitalIntraoralXRayImageStorageForProcessing*  
*CTImageStorage*  
*EnhancedCTImageStorage*  
*UltrasoundMultiframeImageStorageRetired*  
*UltrasoundMultiframeImageStorage*  
*MRIImageStorage*  
*EnhancedMRIImageStorage*  
*MRSpectroscopyStorage*  
*NuclearMedicineImageStorageRetired*  
*UltrasoundImageStorageRetired*  
*UltrasoundImageStorage*  
*SecondaryCaptureImageStorage*  
*MultiframeSingleBitSecondaryCaptureImageStorage*  
*MultiframeGrayscaleByteSecondaryCaptureImageStorage*  
*MultiframeGrayscaleWordSecondaryCaptureImageStorage*  
*MultiframeTrueColorSecondaryCaptureImageStorage*  
*StandaloneOverlayStorageRetired*  
*StandaloneCurveStorageRetired*  
*WaveformStorageTrialRetired*  
*GeneralECGWaveformStorage*  
*AmbulatoryECGWaveformStorage*  
*HemodynamicWaveformStorage*  
*CardiacElectrophysiologyWaveformStorage*

*BasicVoiceAudioWaveformStorage*  
*StandaloneModalityLUTStorageRetired*  
*StandaloneVOILUTStorageRetired*  
*GrayscaleSoftcopyPresentationStateStorageSOPClass*  
*ColorSoftcopyPresentationStateStorageSOPClass*  
*PseudoColorSoftcopyPresentationStateStorageSOPClass*  
*BlendingSoftcopyPresentationStateStorageSOPClass*  
*XRayAngiographicImageStorage*  
*EnhancedXAImageStorage*  
*XRayRadiofluoroscopicImageStorage*  
*EnhancedXRImageStorage*  
*XRay3DAngiographicImageStorage*  
*XRay3DCraniofacialImageStorage*  
*XRayAngiographicBiPlaneImageStorageRetired*  
*NuclearMedicineImageStorage*  
*RawDataStorage*  
*SpatialRegistrationStorage*  
*SpatialFiducialsStorage*  
*DeformableSpatialRegistrationStorage*  
*SegmentationStorage*  
*RealWorldValueMappingStorage*  
*VImageStorageTrialRetired*  
*VMultiframeImageStorageTrialRetired*  
*VLEndoscopicImageStorage*  
*VideoEndoscopicImageStorage*  
*VLMicroscopicImageStorage*  
*VideoMicroscopicImageStorage*  
*VLSlideCoordinatesMicroscopicImageStorage*  
*VLPhotographicImageStorage*  
*VideoPhotographicImageStorage*  
*OphthalmicPhotography8BitImageStorage*  
*OphthalmicPhotography16BitImageStorage*  
*StereometricRelationshipStorage*  
*OphthalmicTomographyImageStorage*  
*TextSRStorageTrialRetired*  
*AudioSRStorageTrialRetired*  
*DetailSRStorageTrialRetired*  
*ComprehensiveSRStorageTrialRetired*  
*BasicTextSRStorage*  
*EnhancedSRStorage*  
*ComprehensiveSRStorage*  
*ProcedureLogStorage*

*MammographyCADSRStorage*  
*KeyObjectSelectionDocumentStorage*  
*ChestCADSRStorage*  
*XRayRadiationDoseSRStorage*  
*EncapsulatedPDFStorage*  
*EncapsulatedCDASStorage*  
*PositronEmissionTomographyImageStorage*  
*StandalonePETCurveStorageRetired*  
*RTImageStorage*  
*RTDoseStorage*  
*RTStructureSetStorage*  
*RTBeamsTreatmentRecordStorage*  
*RTPlanStorage*  
*RTBrachyTreatmentRecordStorage*  
*RTTreatmentSummaryRecordStorage*  
*RTIonPlanStorage*  
*RTIonBeamsTreatmentRecordStorage*  
*PatientRootQueryRetrieveInformationModelFIND*  
*PatientRootQueryRetrieveInformationModelMOVE*  
*PatientRootQueryRetrieveInformationModelGET*  
*StudyRootQueryRetrieveInformationModelFIND*  
*StudyRootQueryRetrieveInformationModelMOVE*  
*StudyRootQueryRetrieveInformationModelGET*  
*PatientStudyOnlyQueryRetrieveInformationModelFINDRetired*  
*PatientStudyOnlyQueryRetrieveInformationModelMOVERetired*  
*PatientStudyOnlyQueryRetrieveInformationModelGETRetired*  
*ModalityWorklistInformationModelFIND*  
*GeneralPurposeWorklistInformationModelFIND*  
*GeneralPurposeScheduledProcedureStepSOPClass*  
*GeneralPurposePerformedProcedureStepSOPClass*  
*GeneralPurposeWorklistManagementMetaSOPClass*  
*InstanceAvailabilityNotificationSOPClass*  
*RTBeamsDeliveryInstructionStorageSupplement74FrozenDraft*  
*RTConventionalMachineVerificationSupplement74FrozenDraft*  
*RTIonMachineVerificationSupplement74FrozenDraft*  
*UnifiedWorklistandProcedureStepServiceClass*  
*UnifiedProcedureStepPushSOPClass*  
*UnifiedProcedureStepWatchSOPClass*  
*UnifiedProcedureStepPullSOPClass*  
*UnifiedProcedureStepEventSOPClass*  
*UnifiedWorklistandProcedureStepSOPInstance*  
*GeneralRelevantPatientInformationQuery*

*BreastImagingRelevantPatientInformationQuery*  
*CardiacRelevantPatientInformationQuery*  
*HangingProtocolStorage*  
*HangingProtocolInformationModelFIND*  
*HangingProtocolInformationModelMOVE*  
*ProductCharacteristicsQuerySOPClass*  
*SubstanceApprovalQuerySOPClass*  
*dicomDeviceName*  
*dicomDescription*  
*dicomManufacturer*  
*dicomManufacturerModelName*  
*dicomSoftwareVersion*  
*dicomVendorData*  
*dicomAETitle*  
*dicomNetworkConnectionReference*  
*dicomApplicationCluster*  
*dicomAssociationInitiator*  
*dicomAssociationAcceptor*  
*dicomHostname*  
*dicomPort*  
*dicomSOPClass*  
*dicomTransferRole*  
*dicomTransferSyntax*  
*dicomPrimaryDeviceType*  
*dicomRelatedDeviceReference*  
*dicomPreferredCalledAETitle*  
*dicomTLSCyphersuite*  
*dicomAuthorizedNodeCertificateReference*  
*dicomThisNodeCertificateReference*  
*dicomInstalled*  
*dicomStationName*  
*dicomDeviceSerialNumber*  
*dicomInstitutionName*  
*dicomInstitutionAddress*  
*dicomInstitutionDepartmentName*  
*dicomIssuerOfPatientID*  
*dicomPreferredCallingAETitle*  
*dicomSupportedCharacterSet*  
*dicomConfigurationRoot*  
*dicomDevicesRoot*  
*dicomUniqueAETitlesRegistryRoot*  
*dicomDevice*

*dicomNetworkAE*  
*dicomNetworkConnection*  
*dicomUniqueAETitle*  
*dicomTransferCapability*  
*VLWholeSlideMicroscopyImageStorage*  
*EnhancedUSVolumeStorage*  
*SurfaceSegmentationStorage*  
*BreastTomosynthesisImageStorage*

### 27.290.3.2 enum gdcmm::UIDs::TSType

Enumerator

*uid\_1\_2\_840\_10008\_1\_1*  
*uid\_1\_2\_840\_10008\_1\_2*  
*uid\_1\_2\_840\_10008\_1\_2\_1*  
*uid\_1\_2\_840\_10008\_1\_2\_1\_99*  
*uid\_1\_2\_840\_10008\_1\_2\_2*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_50*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_51*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_52*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_53*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_54*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_55*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_56*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_57*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_58*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_59*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_60*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_61*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_62*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_63*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_64*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_65*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_66*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_70*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_80*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_81*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_90*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_91*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_92*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_93*

*uid\_1\_2\_840\_10008\_1\_2\_4\_94*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_95*  
*uid\_1\_2\_840\_10008\_1\_2\_4\_100*  
*uid\_1\_2\_840\_10008\_1\_2\_5*  
*uid\_1\_2\_840\_10008\_1\_2\_6\_1*  
*uid\_1\_2\_840\_10008\_1\_2\_6\_2*  
*uid\_1\_2\_840\_10008\_1\_3\_10*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_1*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_2*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_3*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_4*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_5*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_6*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_7*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_8*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_9*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_10*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_11*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_12*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_13*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_14*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_15*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_16*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_17*  
*uid\_1\_2\_840\_10008\_1\_4\_1\_18*  
*uid\_1\_2\_840\_10008\_1\_4\_2\_1*  
*uid\_1\_2\_840\_10008\_1\_4\_2\_2*  
*uid\_1\_2\_840\_10008\_1\_9*  
*uid\_1\_2\_840\_10008\_1\_20\_1*  
*uid\_1\_2\_840\_10008\_1\_20\_1\_1*  
*uid\_1\_2\_840\_10008\_1\_20\_2*  
*uid\_1\_2\_840\_10008\_1\_20\_2\_1*  
*uid\_1\_2\_840\_10008\_1\_40*  
*uid\_1\_2\_840\_10008\_1\_40\_1*  
*uid\_1\_2\_840\_10008\_1\_42*  
*uid\_1\_2\_840\_10008\_1\_42\_1*  
*uid\_1\_2\_840\_10008\_2\_6\_1*  
*uid\_1\_2\_840\_10008\_2\_16\_4*  
*uid\_1\_2\_840\_10008\_3\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_1\_1*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_1\_4*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_2\_1*



*uid\_1\_2\_840\_10008\_3\_1\_2\_3\_1*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_3\_2*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_3\_3*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_3\_4*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_3\_5*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_5\_1*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_5\_4*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_5\_5*  
*uid\_1\_2\_840\_10008\_3\_1\_2\_6\_1*  
*uid\_1\_2\_840\_10008\_4\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_4\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_4\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_9*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_9\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_14*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_15*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_16*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_16\_376*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_17*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_17\_376*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_18*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_18\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_22*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_23*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_24*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_24\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_25*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_26*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_27*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_29*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_30*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_31*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_32*  
*uid\_1\_2\_840\_10008\_5\_1\_1\_33*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2\_1*

*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_5*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_8*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_2\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_3\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_4\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_10*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_20*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_1*

*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_67*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_11*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_22*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_33*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_40*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_50*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_59*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_65*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_67*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_128*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_129*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_5*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_6*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_7*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_8*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_9*

*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_31*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_32\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_32\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_32\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_32*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_33*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_4*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_34\_5*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_37\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_37\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_37\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_38\_1*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_38\_2*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_38\_3*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_41*  
*uid\_1\_2\_840\_10008\_5\_1\_4\_42*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_1*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_2*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_3*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_4*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_5*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_6*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_7*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_8*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_9*  
*uid\_1\_2\_840\_10008\_15\_0\_3\_10*

```

uid_1_2_840_10008_15_0_3_11
uid_1_2_840_10008_15_0_3_12
uid_1_2_840_10008_15_0_3_13
uid_1_2_840_10008_15_0_3_14
uid_1_2_840_10008_15_0_3_15
uid_1_2_840_10008_15_0_3_16
uid_1_2_840_10008_15_0_3_17
uid_1_2_840_10008_15_0_3_18
uid_1_2_840_10008_15_0_3_19
uid_1_2_840_10008_15_0_3_20
uid_1_2_840_10008_15_0_3_21
uid_1_2_840_10008_15_0_3_22
uid_1_2_840_10008_15_0_3_23
uid_1_2_840_10008_15_0_3_24
uid_1_2_840_10008_15_0_3_25
uid_1_2_840_10008_15_0_3_26
uid_1_2_840_10008_15_0_3_27
uid_1_2_840_10008_15_0_3_28
uid_1_2_840_10008_15_0_3_29
uid_1_2_840_10008_15_0_3_30
uid_1_2_840_10008_15_0_3_31
uid_1_2_840_10008_15_0_4_1
uid_1_2_840_10008_15_0_4_2
uid_1_2_840_10008_15_0_4_3
uid_1_2_840_10008_15_0_4_4
uid_1_2_840_10008_15_0_4_5
uid_1_2_840_10008_15_0_4_6
uid_1_2_840_10008_15_0_4_7
uid_1_2_840_10008_15_0_4_8
uid_1_2_840_10008_5_1_4_1_1_77_1_6
uid_1_2_840_10008_5_1_4_1_1_6_2
uid_1_2_840_10008_5_1_4_1_1_66_5
uid_1_2_840_10008_5_1_4_1_1_13_1_3

```

## 27.290.4 Member Function Documentation

### 27.290.4.1 const char\* gdcm::UIDs::GetName ( ) const

When object is Initialize function return the well known name associated with uid return NULL when not initialized

Examples:

[GenerateStandardSOPClasses.cxx](#).

Referenced by `gdcm::operator<<()`.

27.290.4.2 `static unsigned int gdcmm::UIDs::GetNumberOfTransferSyntaxStrings ( ) [static]`

27.290.4.3 `const char* gdcmm::UIDs::GetString ( ) const`

When object is Initialize function return the uid return NULL when not initialized

Examples:

[GenerateStandardSOPClasses.cxx.](#)

Referenced by `gdcmm::operator<<()`.

27.290.4.4 `static const char* const* gdcmm::UIDs::GetTransferSyntaxString ( unsigned int ts ) [static]`

27.290.4.5 `static TransferSyntaxStringsType gdcmm::UIDs::GetTransferSyntaxStrings ( ) [static]`

27.290.4.6 `static const char* gdcmm::UIDs::GetUIDName ( unsigned int ts ) [static]`

27.290.4.7 `static const char* gdcmm::UIDs::GetUIDString ( unsigned int ts ) [static]`

27.290.4.8 `gdcmm::UIDs::operator TSType ( ) const [inline]`

27.290.4.9 `bool gdcmm::UIDs::SetFromUID ( const char * str )`

Initialize object from a string (a uid number) return false on error, and internal state is set to 0

Examples:

[GenerateStandardSOPClasses.cxx.](#)

The documentation for this class was generated from the following file:


- [gdcmmUIDs.h](#)

## 27.291 gdcmm::network::ULAction Class Reference

**ULAction** A [ULConnection](#) in a given ULState can perform certain ULActions. This base class provides the interface for running those ULActions on a given [ULConnection](#).

```
#include <gdcmmULAction.h>
```

Inheritance diagram for gdcmm::network::ULAction:

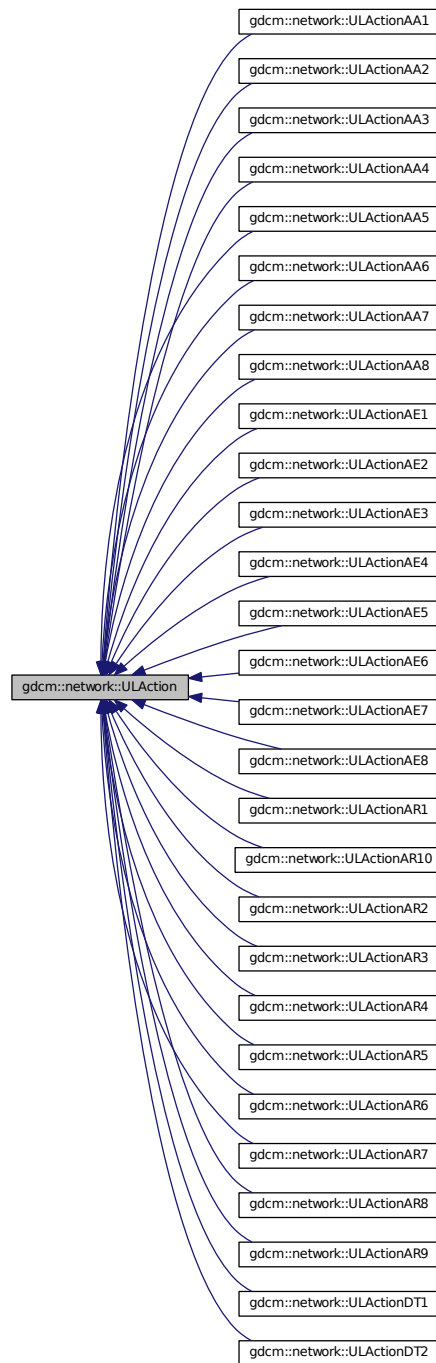

## Public Member Functions

- [ULAction](#) ()
- virtual [~ULAction](#) ()

- virtual [EStateID PerformAction](#) ([Subject \\*s](#), [ULEvent &inEvent](#), [ULConnection &inConnection](#), bool &outWaiting, [ForEvent](#), [EEventID &outRaisedEvent](#))=0

### 27.291.1 Detailed Description

[ULAction](#) A [ULConnection](#) in a given [ULState](#) can perform certain [ULActions](#). This base class provides the interface for running those [ULActions](#) on a given [ULConnection](#).

Essentially, the [ULConnectionManager](#) will take this object, determined from the current [ULState](#) of the [ULConnection](#), and pass the [ULConnection](#) object to the [ULAction](#). The [ULAction](#) will then invoke whatever necessary commands are required by a given action.

The result of a [ULAction](#) is a [ULEvent](#) (ie, what happened as a result of the action).

This [ULEvent](#) is passed to the [ULState](#), so that the transition to the next state can occur.

Actions are associated with Payloads— be those filestreams, AETitles to establish connections, whatever. The actual parameters that the user will pass via an action will come through a Payload object, which should, in itself, be some [gdcm](#)-based object (but not all objects can be payloads; sending a single dataelement as a payload isn't meaningful). As such, each action has its own particular payload.

For the sake of keeping files together, both the particular payload class and the action class will be defined in the same header file. Payloads should JUST be data (or streams), NO METHODS.

Some actions perform changes that should raise events on the local system, and some actions perform changes that will require waiting for events from the remote system.

Therefore, this base action has been modified so that those events are set by each action. When the event loop runs an action, it will then test to see if a local event was raised by the action, and if so, perform the appropriate subsequent action. If the action requires waiting for a response from the remote system, then the event loop will sit there (presumably with the [ARTIM](#) timer running) and wait for a response from the remote system. Once a response is obtained, then the rest of the state transitions can happen.

### 27.291.2 Constructor & Destructor Documentation

27.291.2.1 `gdcm::network::ULAction::ULAction ( ) [inline]`

27.291.2.2 `virtual gdcm::network::ULAction::~~ULAction ( ) [inline],[virtual]`

### 27.291.3 Member Function Documentation

27.291.3.1 `virtual EStateID gdcm::network::ULAction::PerformAction ( Subject \* s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [pure virtual]`

Implemented in [gdcm::network::ULActionAR10](#), [gdcm::network::ULActionAR9](#), [gdcm::network::ULActionAE8](#), [gdcm::network::ULActionAA8](#), [gdcm::network::ULActionAR8](#), [gdcm::network::ULActionAE7](#), [gdcm::network::ULActionAA7](#), [gdcm::network::ULActionAR7](#), [gdcm::network::ULActionAE6](#), [gdcm::network::ULActionAA6](#), [gdcm::network::ULActionAR6](#), [gdcm::network::ULActionAA5](#), [gdcm::network::ULActionAE5](#), [gdcm::network::ULActionAR5](#), [gdcm::network::ULActionAA4](#), [gdcm::network::ULActionAE4](#), [gdcm::network::ULActionAR4](#), [gdcm::network::ULActionAA3](#), [gdcm::network::ULActionAE3](#), [gdcm::network::ULActionAR3](#), [gdcm::network::ULActionAA2](#), [gdcm::network::ULActionAE2](#), [gdcm::network::ULActionAR2](#), [gdcm::network::ULActionDT2](#), [gdcm::network::ULActionAA1](#), [gdcm::network::ULActionAE1](#), [gdcm::network::ULActionAR1](#), and [gdcm::network::ULActionDT1](#).

The documentation for this class was generated from the following file:


- [gdcmULAction.h](#)




## 27.292 gdcmm::network::ULActionAA1 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA1:

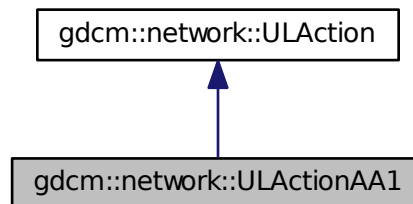

Collaboration diagram for gdcmm::network::ULActionAA1:

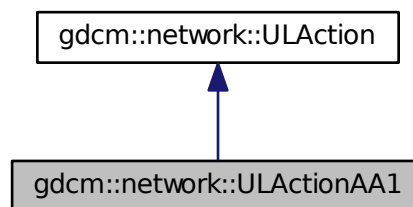

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.292.1 Member Function Documentation

27.292.1.1 `EStateID gdcmm::network::ULActionAA1::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.293 gdcmm::network::ULActionAA2 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA2:

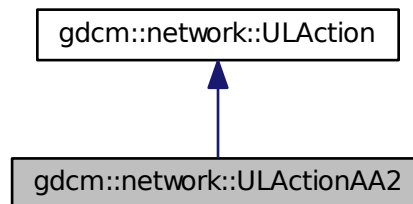

Collaboration diagram for gdcmm::network::ULActionAA2:

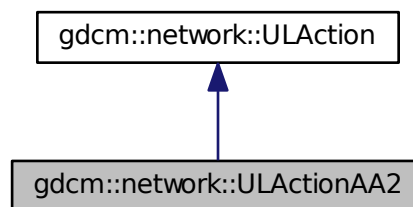

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.293.1 Member Function Documentation

27.293.1.1 [EStateID gdcmm::network::ULActionAA2::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.294 gdcmm::network::ULActionAA3 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA3:

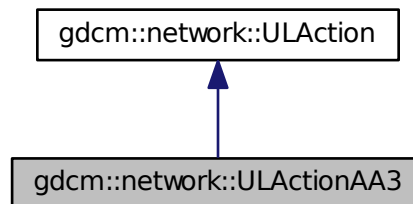

Collaboration diagram for gdcmm::network::ULActionAA3:

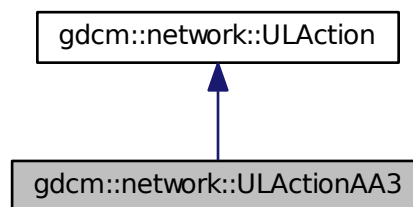

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.294.1 Member Function Documentation

27.294.1.1 `EStateID gdcmm::network::ULActionAA3::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.295 gdcmm::network::ULActionAA4 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA4:

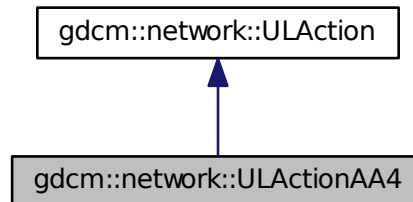

Collaboration diagram for gdcmm::network::ULActionAA4:

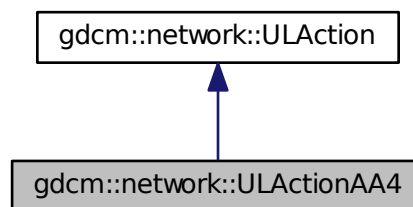

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.295.1 Member Function Documentation

27.295.1.1 [EStateID gdcmm::network::ULActionAA4::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.296 gdcmm::network::ULActionAA5 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA5:

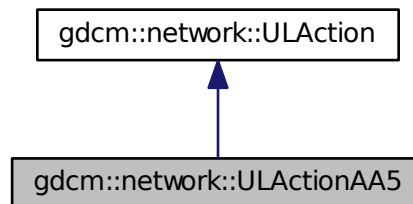

Collaboration diagram for gdcmm::network::ULActionAA5:

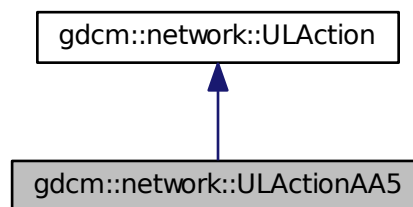

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.296.1 Member Function Documentation

27.296.1.1 [EStateID gdcmm::network::ULActionAA5::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) &  
*inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [virtual]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.297 gdcmm::network::ULActionAA6 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA6:

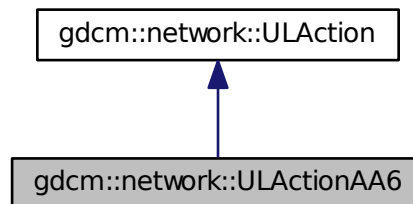

Collaboration diagram for gdcmm::network::ULActionAA6:

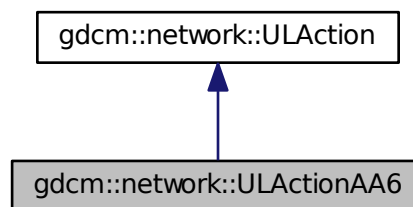

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.297.1 Member Function Documentation

27.297.1.1 `EStateID gdcmm::network::ULActionAA6::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.298 gdcmm::network::ULActionAA7 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA7:

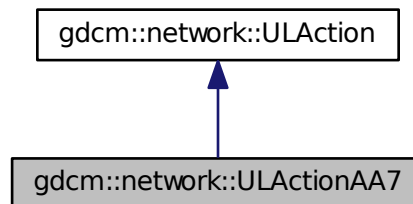

Collaboration diagram for gdcmm::network::ULActionAA7:

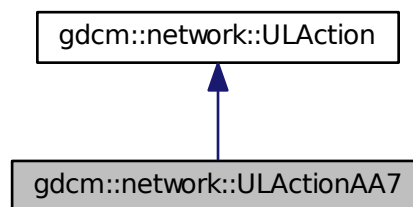

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.298.1 Member Function Documentation

27.298.1.1 [EStateID gdcmm::network::ULActionAA7::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAA.h](#)


## 27.299 gdcmm::network::ULActionAA8 Class Reference

```
#include <gdcmmULActionAA.h>
```

Inheritance diagram for gdcmm::network::ULActionAA8:

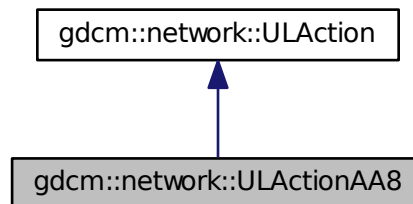

Collaboration diagram for gdcmm::network::ULActionAA8:

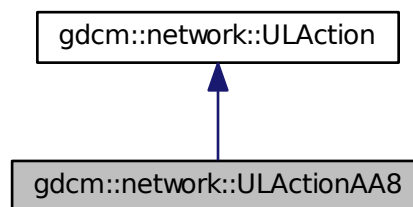

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.299.1 Member Function Documentation

27.299.1.1 [EStateID gdcmm::network::ULActionAA8::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).

The documentation for this class was generated from the following file:


- [gdcmmULActionAA.h](#)




## 27.300 gdcM::network::ULActionAE1 Class Reference

```
#include <gdcMULActionAE.h>
```

Inheritance diagram for gdcM::network::ULActionAE1:

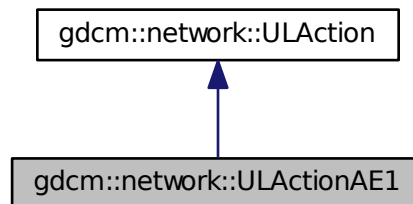

Collaboration diagram for gdcM::network::ULActionAE1:

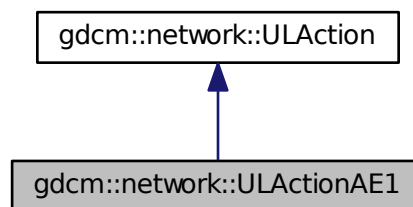

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.300.1 Member Function Documentation

27.300.1.1 [EStateID gdcM::network::ULActionAE1::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) `[virtual]`

Implements [gdcM::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcMULActionAE.h](#)


## 27.301 gdcmm::network::ULActionAE2 Class Reference

```
#include <gdcmmULActionAE.h>
```

Inheritance diagram for gdcmm::network::ULActionAE2:

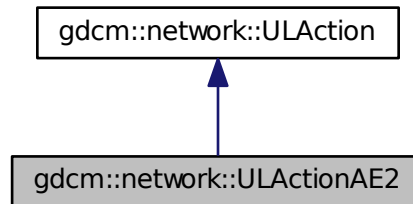

Collaboration diagram for gdcmm::network::ULActionAE2:

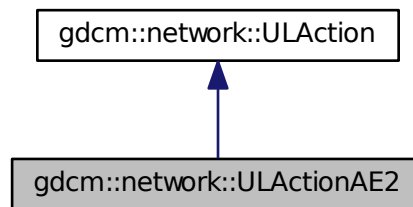

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.301.1 Member Function Documentation

27.301.1.1 [EStateID gdcmm::network::ULActionAE2::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) `[virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAE.h](#)


## 27.302 gdcM::network::ULActionAE3 Class Reference

```
#include <gdcMULActionAE.h>
```

Inheritance diagram for gdcM::network::ULActionAE3:

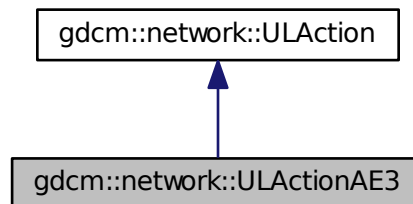

Collaboration diagram for gdcM::network::ULActionAE3:

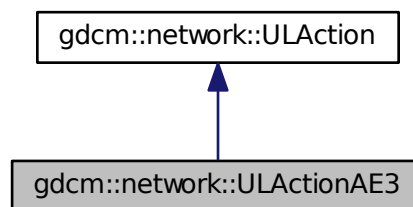

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.302.1 Member Function Documentation

27.302.1.1 [EStateID gdcM::network::ULActionAE3::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [virtual]

Implements [gdcM::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcMULActionAE.h](#)


## 27.303 gdcmm::network::ULActionAE4 Class Reference

```
#include <gdcmmULActionAE.h>
```

Inheritance diagram for gdcmm::network::ULActionAE4:

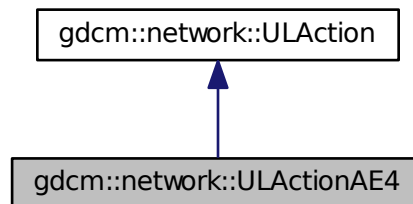

Collaboration diagram for gdcmm::network::ULActionAE4:

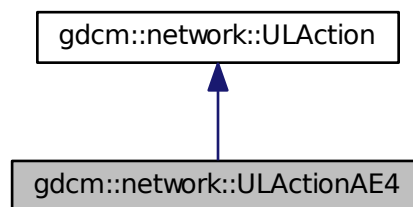

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.303.1 Member Function Documentation

27.303.1.1 `EStateID gdcmm::network::ULActionAE4::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAE.h](#)


## 27.304 gdcM::network::ULActionAE5 Class Reference

```
#include <gdcMULActionAE.h>
```

Inheritance diagram for gdcM::network::ULActionAE5:

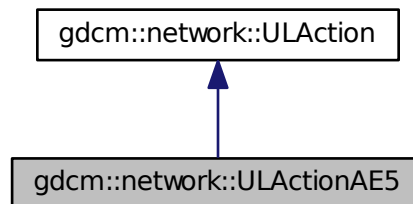

Collaboration diagram for gdcM::network::ULActionAE5:

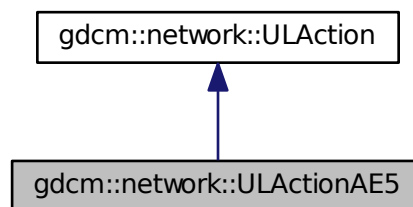

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.304.1 Member Function Documentation

27.304.1.1 [EStateID gdcM::network::ULActionAE5::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) `[virtual]`

Implements [gdcM::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcMULActionAE.h](#)


## 27.305 gdcmm::network::ULActionAE6 Class Reference

```
#include <gdcmmULActionAE.h>
```

Inheritance diagram for gdcmm::network::ULActionAE6:

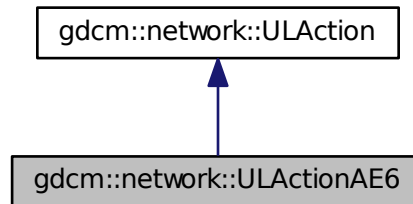

Collaboration diagram for gdcmm::network::ULActionAE6:

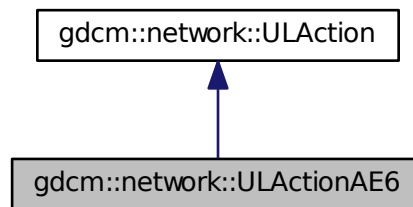

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.305.1 Member Function Documentation

27.305.1.1 `EStateID gdcmm::network::ULActionAE6::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent )` [virtual]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAE.h](#)


## 27.306 gdcmm::network::ULActionAE7 Class Reference

```
#include <gdcmmULActionAE.h>
```

Inheritance diagram for gdcmm::network::ULActionAE7:

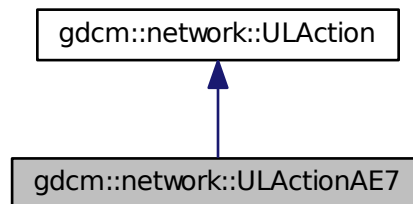

Collaboration diagram for gdcmm::network::ULActionAE7:

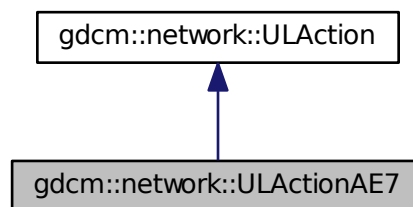

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.306.1 Member Function Documentation

27.306.1.1 [EStateID gdcmm::network::ULActionAE7::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) `[virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAE.h](#)


## 27.307 gdcmm::network::ULActionAE8 Class Reference

```
#include <gdcmmULActionAE.h>
```

Inheritance diagram for gdcmm::network::ULActionAE8:

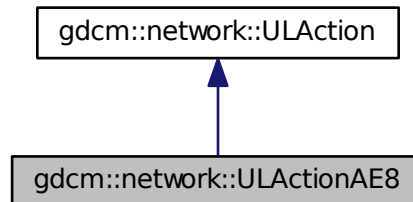

Collaboration diagram for gdcmm::network::ULActionAE8:

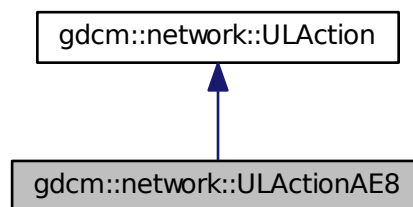

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.307.1 Member Function Documentation

27.307.1.1 [EStateID gdcmm::network::ULActionAE8::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).

The documentation for this class was generated from the following file:


- [gdcmmULActionAE.h](#)




## 27.308 gdcM::network::ULActionAR1 Class Reference

```
#include <gdcMULActionAR.h>
```

Inheritance diagram for gdcM::network::ULActionAR1:

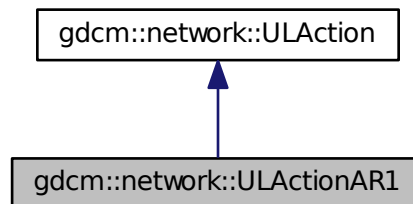

Collaboration diagram for gdcM::network::ULActionAR1:

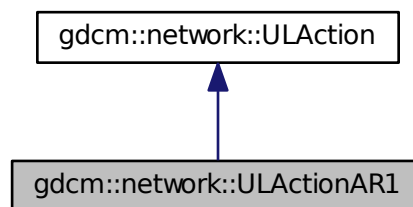

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.308.1 Member Function Documentation

27.308.1.1 [EStateID gdcM::network::ULActionAR1::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [virtual]

Implements [gdcM::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcMULActionAR.h](#)


## 27.309 gdcmm::network::ULActionAR10 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR10:

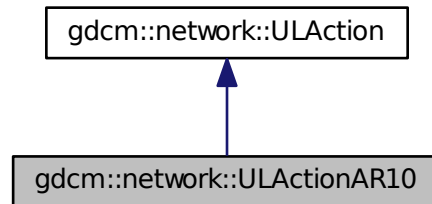

Collaboration diagram for gdcmm::network::ULActionAR10:

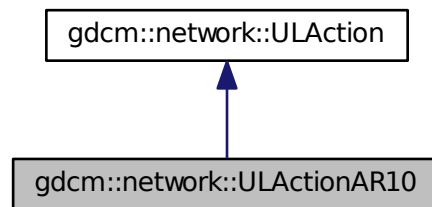

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.309.1 Member Function Documentation

27.309.1.1 **EStateID** `gdcmm::network::ULActionAR10::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent )` [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.310 gdcM::network::ULActionAR2 Class Reference

```
#include <gdcMULActionAR.h>
```

Inheritance diagram for gdcM::network::ULActionAR2:

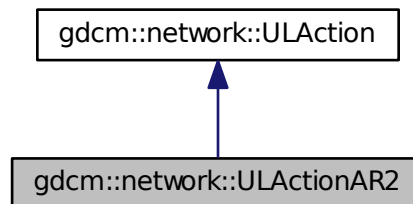

Collaboration diagram for gdcM::network::ULActionAR2:

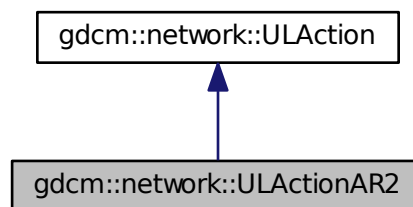

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.310.1 Member Function Documentation

27.310.1.1 [EStateID gdcM::network::ULActionAR2::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [virtual]

Implements [gdcM::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcMULActionAR.h](#)


## 27.311 gdcmm::network::ULActionAR3 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR3:

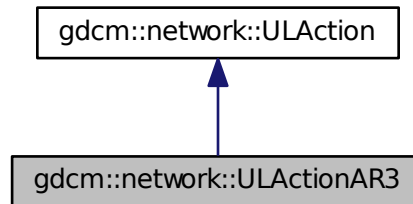

Collaboration diagram for gdcmm::network::ULActionAR3:

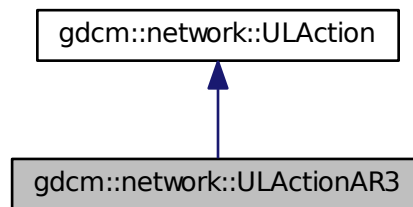

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.311.1 Member Function Documentation

27.311.1.1 [EStateID gdcmm::network::ULActionAR3::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.312 gdcmm::network::ULActionAR4 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR4:

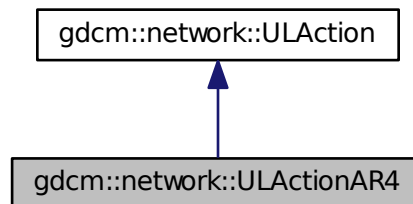

Collaboration diagram for gdcmm::network::ULActionAR4:

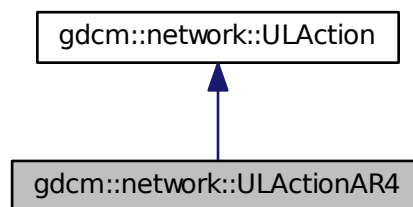

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.312.1 Member Function Documentation

27.312.1.1 `EStateID gdcmm::network::ULActionAR4::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.313 gdcmm::network::ULActionAR5 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR5:

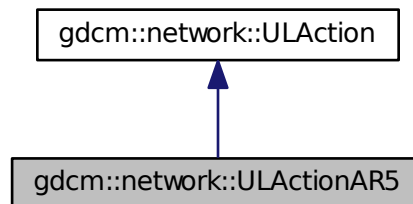

Collaboration diagram for gdcmm::network::ULActionAR5:

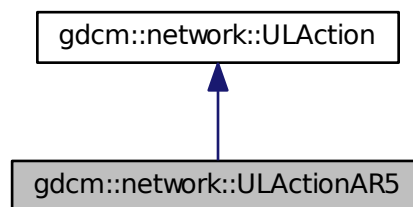

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.313.1 Member Function Documentation

27.313.1.1 `EStateID gdcmm::network::ULActionAR5::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.314 gdcmm::network::ULActionAR6 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR6:

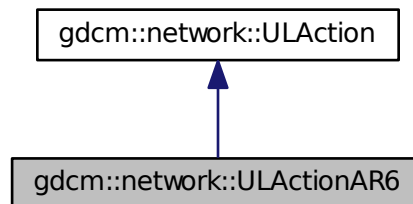

Collaboration diagram for gdcmm::network::ULActionAR6:

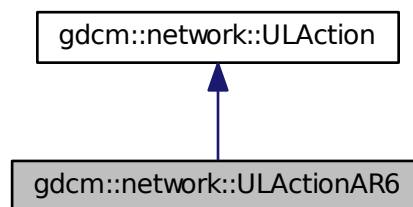

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor↵  
Event, [EEventID](#) &outRaisedEvent)

### 27.314.1 Member Function Documentation

27.314.1.1 [EStateID gdcmm::network::ULActionAR6::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [virtual]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.315 gdcmm::network::ULActionAR7 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR7:

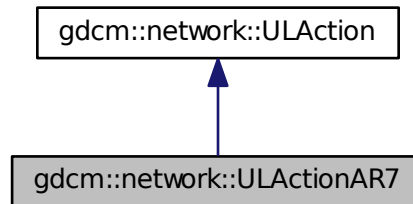

Collaboration diagram for gdcmm::network::ULActionAR7:

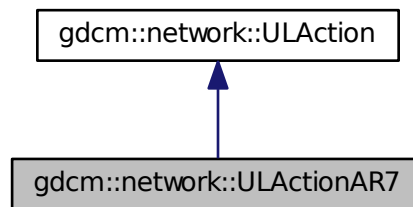

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.315.1 Member Function Documentation

27.315.1.1 `EStateID gdcmm::network::ULActionAR7::PerformAction ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) [virtual]`

Implements [gdcmm::network::ULAction](#).

The documentation for this class was generated from the following file:


- [gdcmmULActionAR.h](#)




## 27.316 gdcmm::network::ULActionAR8 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR8:

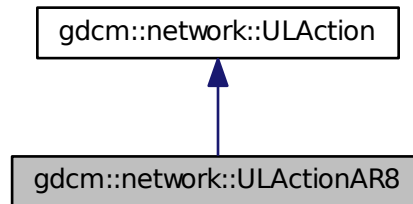

Collaboration diagram for gdcmm::network::ULActionAR8:

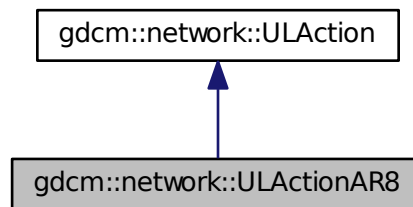

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [EEventID](#) &outRaisedEvent)

### 27.316.1 Member Function Documentation

27.316.1.1 [EStateID gdcmm::network::ULActionAR8::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [\[virtual\]](#)

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.317 gdcmm::network::ULActionAR9 Class Reference

```
#include <gdcmmULActionAR.h>
```

Inheritance diagram for gdcmm::network::ULActionAR9:

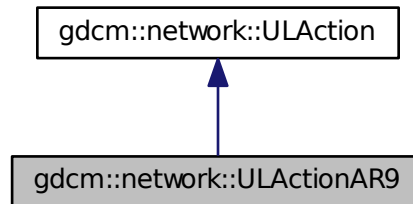

Collaboration diagram for gdcmm::network::ULActionAR9:

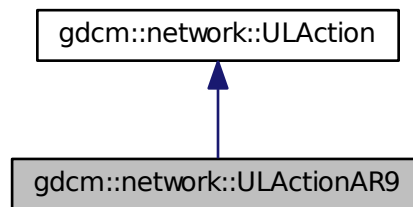

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.317.1 Member Function Documentation

27.317.1.1 [EStateID gdcmm::network::ULActionAR9::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionAR.h](#)


## 27.318 gdcmm::network::ULActionDT1 Class Reference

```
#include <gdcmmULActionDT.h>
```

Inheritance diagram for gdcmm::network::ULActionDT1:

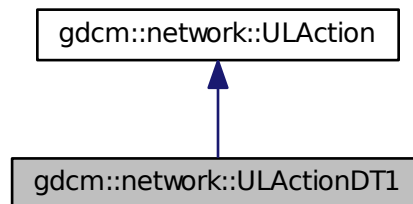

Collaboration diagram for gdcmm::network::ULActionDT1:

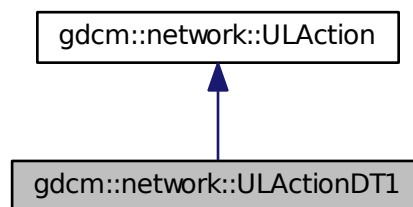

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.318.1 Member Function Documentation

27.318.1.1 [EStateID gdcmm::network::ULActionDT1::PerformAction](#) ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).


The documentation for this class was generated from the following file:

- [gdcmmULActionDT.h](#)


## 27.319 gdcmm::network::ULActionDT2 Class Reference

```
#include <gdcmmULActionDT.h>
```

Inheritance diagram for gdcmm::network::ULActionDT2:

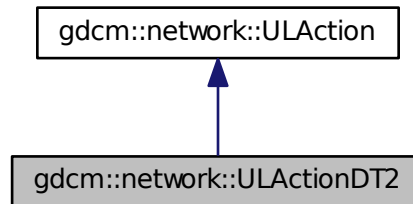

Collaboration diagram for gdcmm::network::ULActionDT2:

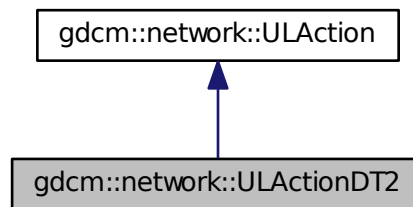

### Public Member Functions

- [EStateID PerformAction](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingFor←  
Event, [EEventID](#) &outRaisedEvent)

### 27.319.1 Member Function Documentation

27.319.1.1 [EStateID](#) `gdcmm::network::ULActionDT2::PerformAction` ( [Subject](#) \* s, [ULEvent](#) & *inEvent*, [ULConnection](#) & *inConnection*, bool & *outWaitingForEvent*, [EEventID](#) & *outRaisedEvent* ) [[virtual](#)]

Implements [gdcmm::network::ULAction](#).

The documentation for this class was generated from the following file:


- [gdcmmULActionDT.h](#)

## 27.320 gdcm::network::ULBasicCallback Class Reference


[ULBasicCallback](#) This is the most basic of callbacks for how the [ULConnectionManager](#) handles incoming datasets. DataSets are just concatenated to the mDataSets vector, and the result can be pulled out of the vector by later code. Alternatives to this method include progress updates, saving to disk, etc. This class is NOT THREAD SAFE. Access the dataset vector after the entire set of datasets has been returned by the [ULConnectionManager](#).

```
#include <gdcmULBasicCallback.h>
```

Inheritance diagram for gdcm::network::ULBasicCallback:

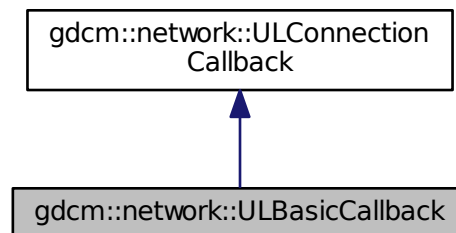

Collaboration diagram for gdcm::network::ULBasicCallback:

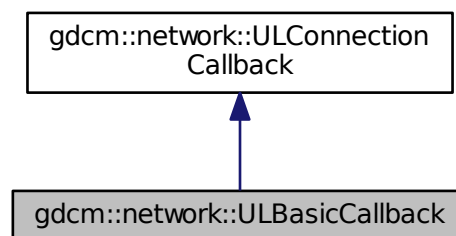

### Public Member Functions

- [ULBasicCallback](#) ()
- virtual [~ULBasicCallback](#) ()
- std::vector< [DataSet](#) > const & [GetDataSets](#) () const
- std::vector< [DataSet](#) > const & [GetResponses](#) () const
- virtual void [HandleDataSet](#) (const [DataSet](#) &inDataSet)
- virtual void [HandleResponse](#) (const [DataSet](#) &inDataSet)

## Additional Inherited Members

### 27.320.1 Detailed Description

[ULBasicCallback](#) This is the most basic of callbacks for how the [ULConnectionManager](#) handles incoming datasets. DataSets are just concatenated to the mDataSets vector, and the result can be pulled out of the vector by later code. Alternatives to this method include progress updates, saving to disk, etc. This class is NOT THREAD SAFE. Access the dataset vector after the entire set of datasets has been returned by the [ULConnectionManager](#).

### 27.320.2 Constructor & Destructor Documentation

27.320.2.1 `gdcm::network::ULBasicCallback::ULBasicCallback ( ) [inline]`

27.320.2.2 `virtual gdcm::network::ULBasicCallback::~~ULBasicCallback ( ) [inline],[virtual]`

### 27.320.3 Member Function Documentation

27.320.3.1 `std::vector<DataSet> const& gdcm::network::ULBasicCallback::GetDataSets ( ) const`

27.320.3.2 `std::vector<DataSet> const& gdcm::network::ULBasicCallback::GetResponses ( ) const`

27.320.3.3 `virtual void gdcm::network::ULBasicCallback::HandleDataSet ( const DataSet & inDataSet ) [virtual]`

Implements [gdcm::network::ULConnectionCallback](#).

27.320.3.4 `virtual void gdcm::network::ULBasicCallback::HandleResponse ( const DataSet & inDataSet ) [virtual]`

Implements [gdcm::network::ULConnectionCallback](#).

The documentation for this class was generated from the following file:

- [gdcmULBasicCallback.h](#)

## 27.321 gdcm::network::ULConnection Class Reference

[ULConnection](#) This is the class that contains the socket to another machine, and passes data through itself, as well as maintaining a sense of state.

```
#include <gdcmULConnection.h>
```

### Public Member Functions

- [ULConnection](#) (const [ULConnectionInfo](#) &inUserInformation)
- virtual [~ULConnection](#) ()
- void [AddAcceptedPresentationContext](#) (const [PresentationContextAC](#) &inPC)
- [PresentationContextRQ FindContext](#) (const [DataElement](#) &de) const
- std::vector< [PresentationContextAC](#) > const & [GetAcceptedPresentationContexts](#) () const
- std::vector< [PresentationContextAC](#) > & [GetAcceptedPresentationContexts](#) ()
- const [ULConnectionInfo](#) & [GetConnectionInfo](#) () const
- uint32\_t [GetMaxPDUSize](#) () const

- const [PresentationContextAC](#) \* [GetPresentationContextACByID](#) (uint8\_t id) const
- uint8\_t [GetPresentationContextIDFromPresentationContext](#) ([PresentationContextRQ](#) const &pc) const  
*return 0 upon error*
- const [PresentationContextRQ](#) \* [GetPresentationContextRQByID](#) (uint8\_t id) const
- std::vector< [PresentationContextRQ](#) > const & [GetPresentationContexts](#) () const
- std::iostream \* [GetProtocol](#) ()
- [EStateID](#) [GetState](#) () const
- [ARTIMTimer](#) & [GetTimer](#) ()
- bool [InitializeConnection](#) ()  
*used to establish scu connections*
- bool [InitializeIncomingConnection](#) ()  
*used to establish scp connections*
- void [SetMaxPDUSize](#) (uint32\_t inSize)
- void [SetPresentationContexts](#) (const std::vector< [PresentationContextRQ](#) > &inContexts)
- void [SetPresentationContexts](#) (const std::vector< [PresentationContext](#) > &inContexts)
- void [SetState](#) (const [EStateID](#) &inState)
- void [StopProtocol](#) ()

## Friends

- class [ULActionAE6](#)
- class [ULConnectionManager](#)

## 27.321.1 Detailed Description

[ULConnection](#) This is the class that contains the socket to another machine, and passes data through itself, as well as maintaining a sense of state.

The [ULConnectionManager](#) tells the [ULConnection](#) what data can actually be sent.

This class is done this way so that it can be eventually be replaced with a [ULSecureConnection](#), if such a protocol is warranted, so that all data that passes through can be managed through a secure connection. For now, this class provides a simple pass-through mechanism to the socket itself.

So, for instance, a gdcm object will be passes to this object, and it will then get passed along the connection, if that connection is in the proper state to do so.

For right now, this class is not directly intended to be inherited from, but the potential for future [ULSecureConnection](#) warrants the addition, rather than having everything be managed from within the [ULConnectionManager](#) (or this class) without a wrapper.

## 27.321.2 Constructor & Destructor Documentation

27.321.2.1 `gdcm::network::ULConnection::ULConnection ( const ULConnectionInfo & inUserInformation )`

27.321.2.2 `virtual gdcm::network::ULConnection::~~ULConnection ( ) [virtual]`

## 27.321.3 Member Function Documentation

27.321.3.1 `void gdcm::network::ULConnection::AddAcceptedPresentationContext ( const PresentationContextAC & inPC )`

27.321.3.2 **PresentationContextRQ** gdcm::network::ULConnection::FindContext ( const DataElement & *de* ) const

27.321.3.3 std::vector<PresentationContextAC> const& gdcm::network::ULConnection::GetAcceptedPresentationContexts ( ) const

27.321.3.4 std::vector<PresentationContextAC>& gdcm::network::ULConnection::GetAcceptedPresentationContexts ( )

27.321.3.5 const ULConnectionInfo& gdcm::network::ULConnection::GetConnectionInfo ( ) const

27.321.3.6 uint32\_t gdcm::network::ULConnection::GetMaxPDUSize ( ) const

27.321.3.7 const PresentationContextAC\* gdcm::network::ULConnection::GetPresentationContextACByID ( uint8\_t *id* ) const

27.321.3.8 uint8\_t gdcm::network::ULConnection::GetPresentationContextIDFromPresentationContext ( PresentationContextRQ const & *pc* ) const

return 0 upon error

27.321.3.9 const PresentationContextRQ\* gdcm::network::ULConnection::GetPresentationContextRQByID ( uint8\_t *id* ) const

27.321.3.10 std::vector<PresentationContextRQ> const& gdcm::network::ULConnection::GetPresentationContexts ( ) const

27.321.3.11 std::iostream\* gdcm::network::ULConnection::GetProtocol ( )

27.321.3.12 EStateID gdcm::network::ULConnection::GetState ( ) const

27.321.3.13 ARTIMTimer& gdcm::network::ULConnection::GetTimer ( )

27.321.3.14 bool gdcm::network::ULConnection::InitializeConnection ( )

used to establish scu connections

27.321.3.15 bool gdcm::network::ULConnection::InitializeIncomingConnection ( )

used to establish scp connections

27.321.3.16 void gdcm::network::ULConnection::SetMaxPDUSize ( uint32\_t *inSize* )

27.321.3.17 void gdcm::network::ULConnection::SetPresentationContexts ( const std::vector< PresentationContextRQ > & *inContexts* )

27.321.3.18 void gdcm::network::ULConnection::SetPresentationContexts ( const std::vector< PresentationContext > & *inContexts* )

27.321.3.19 void gdcm::network::ULConnection::SetState ( const EStateID & *inState* )

27.321.3.20 void gdcm::network::ULConnection::StopProtocol ( )

## 27.321.4 Friends And Related Function Documentation



27.321.4.1 friend class **ULActionAE6** [friend]

27.321.4.2 friend class **ULConnectionManager** [friend]


The documentation for this class was generated from the following file:

- [gdcmULConnection.h](#)

## 27.322 gdcm::network::ULConnectionCallback Class Reference

```
#include <gdcmULConnectionCallback.h>
```

Inheritance diagram for gdcm::network::ULConnectionCallback:

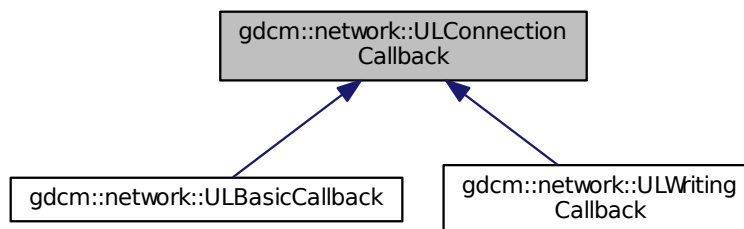

### Public Member Functions

- [ULConnectionCallback](#) ()
- virtual [~ULConnectionCallback](#) ()
- bool [DataSetHandles](#) () const
- virtual void [HandleDataSet](#) (const [DataSet](#) &inDataSet)=0
- virtual void [HandleResponse](#) (const [DataSet](#) &inDataSet)=0
- void [ResetHandledDataSet](#) ()
- void [SetImplicitFlag](#) (const bool imp)

### Protected Member Functions

- void [DataSetHandled](#) ()

### Protected Attributes

- bool [mImplicit](#)

### 27.322.1 Detailed Description

When a dataset comes back from a query/move/etc, the result can either be stored entirely in memory, or could be stored on disk. This class provides a mechanism to indicate what the [ULConnectionManager](#) should do with datasets that are produced through query results. The [ULConnectionManager](#) will call the `HandleDataSet` function during the course of receiving datasets. Particular implementations should fill in what that function does, including updating progress, etc. NOTE: since cmove requires that multiple event loops be employed, the callback function MUST set `mHandledData` Set to true. otherwise, the cmove event loop handler will not know data was received, and proceed to end the loop prematurely.

### 27.322.2 Constructor & Destructor Documentation

27.322.2.1 `gdcm::network::ULConnectionCallback::ULConnectionCallback ( )` `[inline]`

27.322.2.2 `virtual gdcm::network::ULConnectionCallback::~ULConnectionCallback ( )` `[inline]`, `[virtual]`

### 27.322.3 Member Function Documentation

27.322.3.1 `void gdcm::network::ULConnectionCallback::DataSetHandled ( )` `[inline]`, `[protected]`

27.322.3.2 `bool gdcm::network::ULConnectionCallback::DataSetHandles ( )` `const` `[inline]`

27.322.3.3 `virtual void gdcm::network::ULConnectionCallback::HandleDataSet ( const DataSet & inDataSet )` `[pure virtual]`

Implemented in [gdcm::network::ULBasicCallback](#), and [gdcm::network::ULWritingCallback](#).

27.322.3.4 `virtual void gdcm::network::ULConnectionCallback::HandleResponse ( const DataSet & inDataSet )` `[pure virtual]`

Implemented in [gdcm::network::ULBasicCallback](#), and [gdcm::network::ULWritingCallback](#).

27.322.3.5 `void gdcm::network::ULConnectionCallback::ResetHandledDataSet ( )` `[inline]`

27.322.3.6 `void gdcm::network::ULConnectionCallback::SetImplicitFlag ( const bool imp )` `[inline]`

### 27.322.4 Member Data Documentation

27.322.4.1 `bool gdcm::network::ULConnectionCallback::mImplicit` `[protected]`

The documentation for this class was generated from the following file:

- [gdcmULConnectionCallback.h](#)

## 27.323 gdcm::network::ULConnectionInfo Class Reference

[ULConnectionInfo](#) this class contains all the information about a particular connection as established by the user. That is, it's: User Information Calling AE Title Called AE Title IP address/computer name IP Port A connection must be established with this information, that's subsequently placed into various primitives for actual communication.

```
#include <gdcmULConnectionInfo.h>
```

## Public Member Functions

- [ULConnectionInfo](#) ()
- const char \* [GetCalledAETitle](#) () const
- std::string [GetCalledComputerName](#) () const
- unsigned long [GetCalledIPAddress](#) () const
- int [GetCalledIPPort](#) () const
- const char \* [GetCallingAETitle](#) () const
- unsigned long [GetMaxPDULength](#) () const
- bool [Initialize](#) ([UserInfo](#) const &inUserInfo, const char \*inCalledAETitle, const char \*inCallingAETitle, unsigned long inCalledIPAddress, int inCalledIPPort, std::string inCalledComputerName)
- void [SetMaxPDULength](#) (unsigned long inMaxPDULength)

### 27.323.1 Detailed Description

[ULConnectionInfo](#) this class contains all the information about a particular connection as established by the user. That is, it's: User Information Calling AE Title Called AE Title IP address/computer name IP Port A connection must be established with this information, that's subsequently placed into various primitives for actual communication.

### 27.323.2 Constructor & Destructor Documentation

27.323.2.1 `gdcm::network::ULConnectionInfo::ULConnectionInfo ( )`

### 27.323.3 Member Function Documentation

27.323.3.1 `const char* gdcm::network::ULConnectionInfo::GetCalledAETitle ( ) const`

27.323.3.2 `std::string gdcm::network::ULConnectionInfo::GetCalledComputerName ( ) const`

27.323.3.3 `unsigned long gdcm::network::ULConnectionInfo::GetCalledIPAddress ( ) const`

27.323.3.4 `int gdcm::network::ULConnectionInfo::GetCalledIPPort ( ) const`

27.323.3.5 `const char* gdcm::network::ULConnectionInfo::GetCallingAETitle ( ) const`

27.323.3.6 `unsigned long gdcm::network::ULConnectionInfo::GetMaxPDULength ( ) const`

27.323.3.7 `bool gdcm::network::ULConnectionInfo::Initialize ( UserInfo const & inUserInfo, const char * inCalledAETitle, const char * inCallingAETitle, unsigned long inCalledIPAddress, int inCalledIPPort, std::string inCalledComputerName )`

27.323.3.8 `void gdcm::network::ULConnectionInfo::SetMaxPDULength ( unsigned long inMaxPDULength )`

The documentation for this class was generated from the following file:


- [gdcmULConnectionInfo.h](#)

## 27.324 gdcm::network::ULConnectionManager Class Reference


[ULConnectionManager](#) The [ULConnectionManager](#) performs actions on the [ULConnection](#) given inputs from the user and from the state of what's going on around the connection (ie, timeouts of the ARTIM timer, responses from the peer across the connection, etc).

```
#include <gdcmULConnectionManager.h>
```

Inheritance diagram for `gdcm::network::ULConnectionManager`:

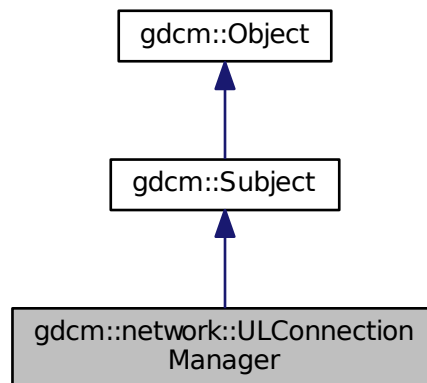

Collaboration diagram for `gdcm::network::ULConnectionManager`:

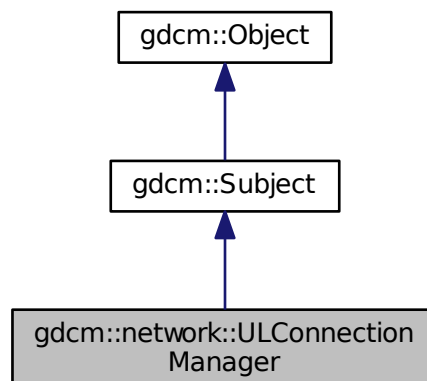

## Public Member Functions

- [ULConnectionManager](#) ()
- [~ULConnectionManager](#) ()
- bool [BreakConnection](#) (const double &inTimeout)
- void [BreakConnectionNow](#) ()
- bool [EstablishConnection](#) (const std::string &inAETitle, const std::string &inConnectAETitle, const std::string &inComputerName, long inIPAddress, uint16\_t inConnectPort, double inTimeout, std::vector< [PresentationContext](#) > const &pcVector)
- bool [EstablishConnectionMove](#) (const std::string &inAETitle, const std::string &inConnectAETitle, const std::string &inComputerName, long inIPAddress, uint16\_t inConnectPort, double inTimeout, uint16\_t inReturnPort, std::vector< [PresentationContext](#) > const &pcVector)
- std::vector< [PresentationDataValue](#) > [SendEcho](#) ()
- std::vector< [DataSet](#) > [SendFind](#) (const [BaseRootQuery](#) \*inRootQuery)
- void [SendFind](#) (const [BaseRootQuery](#) \*inRootQuery, [ULConnectionCallback](#) \*inCallback)
- std::vector< [DataSet](#) > [SendMove](#) (const [BaseRootQuery](#) \*inRootQuery)
- bool [SendMove](#) (const [BaseRootQuery](#) \*inRootQuery, [ULConnectionCallback](#) \*inCallback)  
*return false upon error*
- std::vector< [DataSet](#) > [SendStore](#) (const [File](#) &file)
- void [SendStore](#) (const [File](#) &file, [ULConnectionCallback](#) \*inCallback)  
*callback based API*

## Additional Inherited Members

### 27.324.1 Detailed Description

[ULConnectionManager](#) The [ULConnectionManager](#) performs actions on the [ULConnection](#) given inputs from the user and from the state of what's going on around the connection (ie, timeouts of the ARTIM timer, responses from the peer across the connection, etc).

Its inputs are ULEvents, and it performs ULActions.

### 27.324.2 Constructor & Destructor Documentation

27.324.2.1 `gdcmm::network::ULConnectionManager::ULConnectionManager ( )`

27.324.2.2 `gdcmm::network::ULConnectionManager::~~ULConnectionManager ( )`

### 27.324.3 Member Function Documentation

27.324.3.1 `bool gdcmm::network::ULConnectionManager::BreakConnection ( const double & inTimeout )`

27.324.3.2 `void gdcmm::network::ULConnectionManager::BreakConnectionNow ( )`

27.324.3.3 `bool gdcmm::network::ULConnectionManager::EstablishConnection ( const std::string & inAETitle, const std::string & inConnectAETitle, const std::string & inComputerName, long inIPAddress, uint16_t inConnectPort, double inTimeout, std::vector< PresentationContext > const & pcVector )`

returns true if a connection of the given AETitle (ie, 'this' program) is able to connect to the given AETitle and Port in a certain amount of time providing the connection type will establish the proper exchange syntax with a server; if a different functionality is required, a different connection should be established. returns false if the connection type is 'move'— have to give a return port for move to work as specified.

27.324.3.4 `bool gdcmm::network::ULConnectionManager::EstablishConnectionMove ( const std::string & inAETitle, const std::string & inConnectAETitle, const std::string & inComputerName, long inIPAddress, uint16_t inConnectPort, double inTimeout, uint16_t inReturnPort, std::vector< PresentationContext > const & pcVector )`

returns true for above reasons, but contains the special 'move' port

27.324.3.5 `std::vector<PresentationDataValue> gdcmm::network::ULConnectionManager::SendEcho ( )`

27.324.3.6 `std::vector<DataSet> gdcmm::network::ULConnectionManager::SendFind ( const BaseRootQuery * inRootQuery )`

27.324.3.7 `void gdcmm::network::ULConnectionManager::SendFind ( const BaseRootQuery * inRootQuery, ULConnectionCallback * inCallback )`

27.324.3.8 `std::vector<DataSet> gdcmm::network::ULConnectionManager::SendMove ( const BaseRootQuery * inRootQuery )`

27.324.3.9 `bool gdcmm::network::ULConnectionManager::SendMove ( const BaseRootQuery * inRootQuery, ULConnectionCallback * inCallback )`

return false upon error

27.324.3.10 `std::vector<DataSet> gdcmm::network::ULConnectionManager::SendStore ( const File & file )`

27.324.3.11 `void gdcmm::network::ULConnectionManager::SendStore ( const File & file, ULConnectionCallback * inCallback )`

callback based API

The documentation for this class was generated from the following file:

- [gdcmmULConnectionManager.h](#)

## 27.325 gdcmm::network::ULError Class Reference

[ULError](#) base class for network events.

```
#include <gdcmmULError.h>
```

### Public Member Functions

- [ULError](#) (const [EEventID](#) &inEventID, std::vector< [BasePDU](#) \* > const &inBasePDU)
- [ULError](#) (const [EEventID](#) &inEventID, [BasePDU](#) \*inBasePDU)
- [~ULError](#) ()
- [EEventID](#) [GetEvent](#) () const
- std::vector< [BasePDU](#) \* > const & [GetPDUs](#) () const
- void [SetEvent](#) (const [EEventID](#) &inEvent)
- void [SetPDU](#) (std::vector< [BasePDU](#) \* > const &inPDU)

### 27.325.1 Detailed Description

[ULError](#) base class for network events.

An event consists of the event ID and the data associated with that event.

Note that once a PDU is created, it is now the responsibility of the associated event to destroy it!

### 27.325.2 Constructor & Destructor Documentation

27.325.2.1 `gdcmm::network::ULEvent::ULEvent ( const EEventID & inEventID, std::vector< BasePDU * > const & inBasePDU )` `[inline]`

27.325.2.2 `gdcmm::network::ULEvent::ULEvent ( const EEventID & inEventID, BasePDU * inBasePDU )` `[inline]`

27.325.2.3 `gdcmm::network::ULEvent::~~ULEvent ( )` `[inline]`

### 27.325.3 Member Function Documentation

27.325.3.1 `EEventID gdcmm::network::ULEvent::GetEvent ( ) const` `[inline]`

27.325.3.2 `std::vector<BasePDU*> const& gdcmm::network::ULEvent::GetPDUs ( ) const` `[inline]`

27.325.3.3 `void gdcmm::network::ULEvent::SetEvent ( const EEventID & inEvent )` `[inline]`

27.325.3.4 `void gdcmm::network::ULEvent::SetPDU ( std::vector< BasePDU * > const & inPDU )` `[inline]`

The documentation for this class was generated from the following file:

- [gdcmmULEvent.h](#)

## 27.326 gdcmm::network::ULTransitionTable Class Reference

[ULTransitionTable](#) The transition table of all the ULEvents, new ULActions, and ULStates.

```
#include <gdcmmULTransitionTable.h>
```

### Public Member Functions

- [ULTransitionTable](#) ()
- void [HandleEvent](#) ([Subject](#) \*s, [ULEvent](#) &inEvent, [ULConnection](#) &inConnection, bool &outWaitingForEvent, [E↔EventID](#) &outRaisedEvent) const
- void [PrintTable](#) () const

### 27.326.1 Detailed Description

[ULTransitionTable](#) The transition table of all the ULEvents, new ULActions, and ULStates.

Based roughly on the solutions in `player2.cpp` in the boost examples and this so question: <http://stackoverflow.com/questions/1647631/c-state-machine-design>

The transition table is constructed of `TableRow`s. Each row is based on an event, and an event handler in the `Transition↔Table` object takes a given event, and then finds the given row.

Then, given the current state of the connection, determines the appropriate action to take and then the state to transition to next.

## 27.326.2 Constructor & Destructor Documentation

27.326.2.1 `gdcm::network::ULTransitionTable::ULTransitionTable ( )`

## 27.326.3 Member Function Documentation

27.326.3.1 `void gdcm::network::ULTransitionTable::HandleEvent ( Subject * s, ULEvent & inEvent, ULConnection & inConnection, bool & outWaitingForEvent, EEventID & outRaisedEvent ) const`

27.326.3.2 `void gdcm::network::ULTransitionTable::PrintTable ( ) const`


The documentation for this class was generated from the following file:

- [gdcmULTransitionTable.h](#)

## 27.327 gdcm::network::ULWritingCallback Class Reference


```
#include <gdcmULWritingCallback.h>
```

Inheritance diagram for `gdcm::network::ULWritingCallback`:

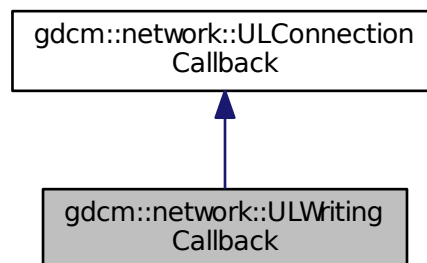



Collaboration diagram for gdcm::network::ULWritingCallback:

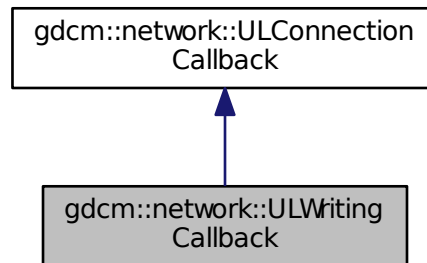

## Public Member Functions

- [ULWritingCallback](#) ()
- virtual [~ULWritingCallback](#) ()
- virtual void [HandleDataSet](#) (const [DataSet](#) &inDataSet)
- virtual void [HandleResponse](#) (const [DataSet](#) &inDataSet)
- void [SetDirectory](#) (const std::string &inDirectoryName)

*provide the directory into which all files are written.*

## Additional Inherited Members

### 27.327.1 Constructor & Destructor Documentation

27.327.1.1 `gdcm::network::ULWritingCallback::ULWritingCallback ( )` `[inline]`

27.327.1.2 `virtual gdcm::network::ULWritingCallback::~~ULWritingCallback ( )` `[inline]`, `[virtual]`

### 27.327.2 Member Function Documentation

27.327.2.1 `virtual void gdcm::network::ULWritingCallback::HandleDataSet ( const DataSet & inDataSet )` `[virtual]`

Implements [gdcm::network::ULConnectionCallback](#).

27.327.2.2 `virtual void gdcm::network::ULWritingCallback::HandleResponse ( const DataSet & inDataSet )` `[virtual]`

Implements [gdcm::network::ULConnectionCallback](#).

27.327.2.3 `void gdcm::network::ULWritingCallback::SetDirectory ( const std::string & inDirectoryName )` `[inline]`

*provide the directory into which all files are written.*

The documentation for this class was generated from the following file:


- [gdcmlWritingCallback.h](#)

## 27.328 gdcml::UNExplicitDataElement Class Reference


Class to read/write a [DataElement](#) as UNExplicit Data [Element](#).

```
#include <gdcmlUNExplicitDataElement.h>
```

Inheritance diagram for gdcml::UNExplicitDataElement:

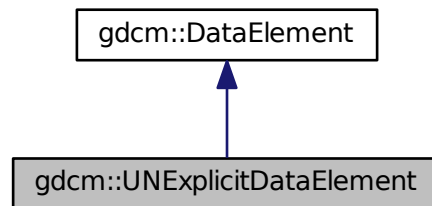

Collaboration diagram for gdcml::UNExplicitDataElement:

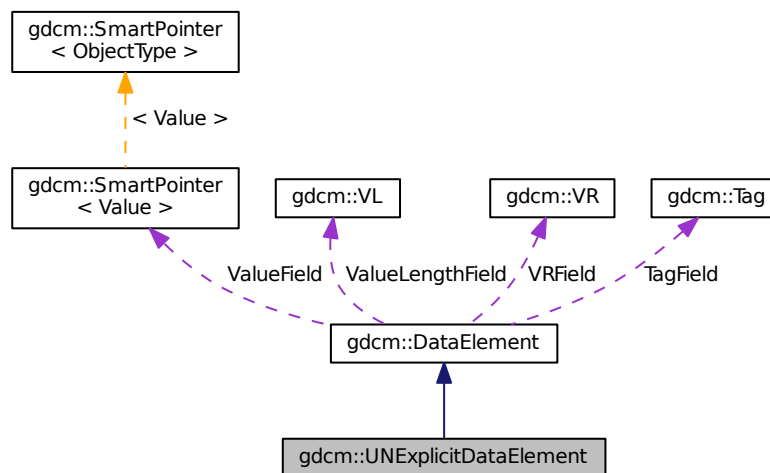

### Public Member Functions

- [VL GetLength](#) () const

- `template<typename TSwap >`  
`std::istream & Read (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & ReadPreValue (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & ReadValue (std::istream &is, bool readvalues=true)`
- `template<typename TSwap >`  
`std::istream & ReadWithLength (std::istream &is, VL &length)`

## Additional Inherited Members

### 27.328.1 Detailed Description

Class to read/write a [DataElement](#) as UNExplicit Data [Element](#).

#### Note

bla

### 27.328.2 Member Function Documentation

27.328.2.1 `VL gdcm::UNExplicitDataElement::GetLength ( ) const`

27.328.2.2 `template<typename TSwap > std::istream& gdcm::UNExplicitDataElement::Read ( std::istream & is )`

27.328.2.3 `template<typename TSwap > std::istream& gdcm::UNExplicitDataElement::ReadPreValue ( std::istream & is )`

27.328.2.4 `template<typename TSwap > std::istream& gdcm::UNExplicitDataElement::ReadValue ( std::istream & is, bool readvalues = true )`

27.328.2.5 `template<typename TSwap > std::istream& gdcm::UNExplicitDataElement::ReadWithLength ( std::istream & is, VL & length )`

The documentation for this class was generated from the following file:


- [gdcmUNExplicitDataElement.h](#)

## 27.329 gdcm::UNExplicitImplicitDataElement Class Reference


Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#) This class gather two known bugs:

```
#include <gdcmUNExplicitImplicitDataElement.h>
```

Inheritance diagram for `gdc::UNExplicitImplicitDataElement`:

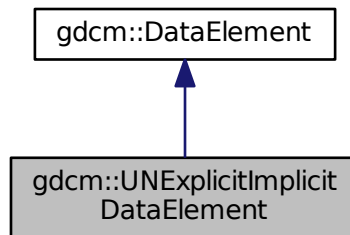

Collaboration diagram for `gdc::UNExplicitImplicitDataElement`:

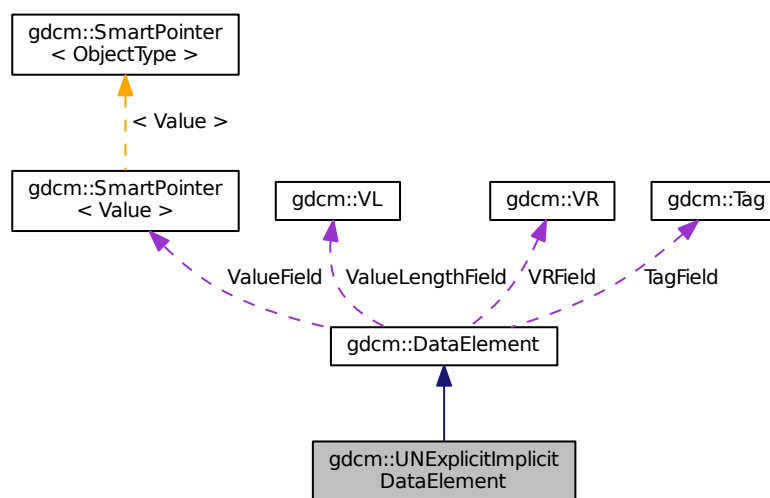

## Public Member Functions

- [VL GetLength](#) () const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadPreValue](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadValue](#) (std::istream &is)

## Additional Inherited Members

### 27.329.1 Detailed Description

Class to read/write a [DataElement](#) as ExplicitImplicit Data [Element](#) This class gather two known bugs:

1. GDCM 1.2.0 would rewrite [VR](#)=UN [Value](#) Length on 2 bytes instead of 4 bytes
2. GDCM 1.2.0 would also rewrite [DataElement](#) as Implicit when the [VR](#) would not be known this would only happen in some very rare cases. gdcm 2.X design could handle bug #1 or #2 exclusively, this class can now handle file which have both issues. See: gdcmData/TherapysGDCM120Bug.dcm

### 27.329.2 Member Function Documentation

27.329.2.1 VL gdcm::UNExplicitImplicitDataElement::GetLength ( ) const

27.329.2.2 template<typename TSwap > std::istream& gdcm::UNExplicitImplicitDataElement::Read ( std::istream & is )

27.329.2.3 template<typename TSwap > std::istream& gdcm::UNExplicitImplicitDataElement::ReadPreValue ( std::istream & is )

27.329.2.4 template<typename TSwap > std::istream& gdcm::UNExplicitImplicitDataElement::ReadValue ( std::istream & is )

The documentation for this class was generated from the following file:

- [gdcmUNExplicitImplicitDataElement.h](#)

## 27.330 gdcm::Unpacker12Bits Class Reference

Pack/Unpack 12 bits pixel into 16bits.

```
#include <gdcmUnpacker12Bits.h>
```

### Static Public Member Functions

- static bool [Pack](#) (char \*out, const char \*in, size\_t n)
- static bool [Unpack](#) (char \*out, const char \*in, size\_t n)

### 27.330.1 Detailed Description

Pack/Unpack 12 bits pixel into 16bits.

- You can only pack an even number of 16bits, which means a multiple of 4 (expressed in bytes)
- You can only unpack a multiple of 3 bytes

This class has no purpose in general purpose DICOM implementation. However to be able to cope with some early ACR-NEMA file generated by a well-known private vendor, one would need to unpack 12bits Stored Pixel [Value](#) into a more standard 16bits Stored Pixel [Value](#).

See also

[Rescaler](#)

## 27.330.2 Member Function Documentation

**27.330.2.1** `static bool gdcm::Unpacker12Bits::Pack ( char * out, const char * in, size_t n )` `[static]`

Pack an array of 16bits where all values are 12bits into a pack form. *n* is the length in bytes of array *in*, *out* will be a fake 8bits array of size  $(n / 2) * 3$

**27.330.2.2** `static bool gdcm::Unpacker12Bits::Unpack ( char * out, const char * in, size_t n )` `[static]`

Unpack an array of 'packed' 12bits data into a more conventional 16bits array. *n* is the length in bytes of array *in*, *out* will be a 16bits array of size  $(n / 3) * 2$

The documentation for this class was generated from the following file:

- [gdcmUnpacker12Bits.h](#)

## 27.331 gdcm::Usage Class Reference

[Usage.](#)

```
#include <gdcmUsage.h>
```

### Public Types

- enum [UsageType](#) {  
[Mandatory](#),  
[Conditional](#),  
[UserOption](#),  
[Invalid](#) }

### Public Member Functions

- [Usage](#) ([UsageType](#) type=[Invalid](#))
- [operator UsageType](#) () const

### Static Public Member Functions

- static const char \* [GetUsageString](#) ([UsageType](#) type)
- static [UsageType](#) [GetUsageType](#) (const char \*type)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [Usage](#) &vr)

### 27.331.1 Detailed Description

#### Usage.

#### Note

A.1.3 IOD Module Table and Functional Group Macro Table This Section of each IOD defines in a tabular form the Modules comprising the IOD. The following information must be specified for each Module in the table:

- The name of the Module or Functional Group
  - A reference to the Section in Annex C which defines the Module or Functional Group
  - The usage of the Module or Functional Group; whether it is:
    - Mandatory (see A.1.3.1) , abbreviated M
    - Conditional (see A.1.3.2) , abbreviated C
    - User Option (see A.1.3.3) , abbreviated U
- The Modules referenced are defined in Annex C. A.1.3.1 MANDATORY MODULES For each IOD, Mandatory Modules shall be supported per the definitions, semantics and requirements defined in Annex C.

A.1.3.2 CONDITIONAL MODULES Conditional Modules are Mandatory Modules if specific conditions are met. If the specified conditions are not met, this Module shall not be supported; that is, no information defined in that Module shall be sent. A.1.3.3 USER OPTION MODULES User Option Modules may or may not be supported. If an optional Module is supported, the Attribute Types specified in the Modules in Annex C shall be supported.

### 27.331.2 Member Enumeration Documentation

#### 27.331.2.1 enum gdcmm::Usage::UsageType

#### Enumerator

***Mandatory***

***Conditional***

***UserOption***

***Invalid***

### 27.331.3 Constructor & Destructor Documentation

27.331.3.1 gdcmm::Usage ( UsageType type = Invalid ) [inline]

### 27.331.4 Member Function Documentation

27.331.4.1 static const char\* gdcmm::Usage::GetUsageString ( UsageType type ) [static]

Referenced by gdcmm::operator<<().

27.331.4.2 static UsageType gdcmm::Usage::GetUsageType ( const char \* type ) [static]

27.331.4.3 gdcmm::Usage::operator UsageType ( ) const [inline]

### 27.331.5 Friends And Related Function Documentation

27.331.5.1 `std::ostream& operator<< ( std::ostream & os, const Usage & vr )` [*friend*]


The documentation for this class was generated from the following file:

- [gdcmUsage.h](#)


## 27.332 gdcm::UserEvent Class Reference

```
#include <gdcmEvent.h>
```

Inheritance diagram for `gdcm::UserEvent`:

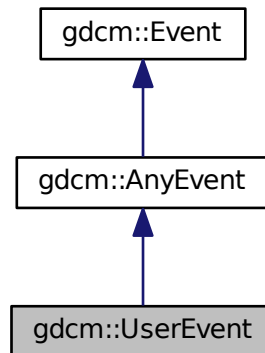

Collaboration diagram for `gdcm::UserEvent`:

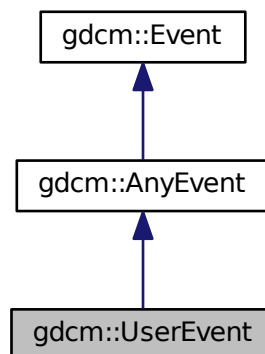



## Additional Inherited Members

The documentation for this class was generated from the following file:

- [gdcmmEvent.h](#)

## 27.333 gdcmm::network::UserInformation Class Reference

[UserInformation Table](#) 9-16 USER INFORMATION ITEM FIELDS.

```
#include <gdcmmUserInformation.h>
```

### Public Member Functions

- [UserInformation](#) ()
- [~UserInformation](#) ()
- void [AddRoleSelectionSub](#) ([RoleSelectionSub](#) const &r)
- void [AddSOPClassExtendedNegociationSub](#) ([SOPClassExtendedNegociationSub](#) const &s)
- const [MaximumLengthSub](#) & [GetMaximumLengthSub](#) () const
- [MaximumLengthSub](#) & [GetMaximumLengthSub](#) ()
- [UserInformation](#) & [operator=](#) (const [UserInformation](#) &)
- void [Print](#) (std::ostream &os) const
- std::istream & [Read](#) (std::istream &is)
- size\_t [Size](#) () const
- const std::ostream & [Write](#) (std::ostream &os) const

### 27.333.1 Detailed Description

[UserInformation Table](#) 9-16 USER INFORMATION ITEM FIELDS.

TODO what is the goal of :

[Table](#) 9-20 USER INFORMATION ITEM FIELDS

### 27.333.2 Constructor & Destructor Documentation

27.333.2.1 [gdcmm::network::UserInformation::UserInformation](#) ( )

27.333.2.2 [gdcmm::network::UserInformation::~~UserInformation](#) ( )

### 27.333.3 Member Function Documentation

27.333.3.1 void [gdcmm::network::UserInformation::AddRoleSelectionSub](#) ( [RoleSelectionSub](#) const & r )

27.333.3.2 void [gdcmm::network::UserInformation::AddSOPClassExtendedNegociationSub](#) ( [SOPClassExtendedNegociationSub](#) const & s )

27.333.3.3 const [MaximumLengthSub](#)& [gdcmm::network::UserInformation::GetMaximumLengthSub](#) ( ) const [inline]

27.333.3.4 **MaxLengthSub**& gdcmm::network::UserInformation::GetMaxLengthSub ( ) [inline]

27.333.3.5 **UserInformation**& gdcmm::network::UserInformation::operator= ( const UserInformation & )

27.333.3.6 void gdcmm::network::UserInformation::Print ( std::ostream & os ) const

27.333.3.7 std::istream& gdcmm::network::UserInformation::Read ( std::istream & is )

27.333.3.8 size\_t gdcmm::network::UserInformation::Size ( ) const

27.333.3.9 const std::ostream& gdcmm::network::UserInformation::Write ( std::ostream & os ) const

The documentation for this class was generated from the following file:

- [gdcmmUserInformation.h](#)

## 27.334 gdcmm::UUIDGenerator Class Reference

Class for generating unique UUID generate DCE 1.1 uid.

```
#include <gdcmmUUIDGenerator.h>
```

### Public Member Functions

- const char \* [Generate](#) ()

### Static Public Member Functions

- static bool [IsValid](#) (const char \*uid)  
*Find out if the string is a valid UUID or not.*

#### 27.334.1 Detailed Description

Class for generating unique UUID generate DCE 1.1 uid.

#### 27.334.2 Member Function Documentation

27.334.2.1 const char\* gdcmm::UUIDGenerator::Generate ( )

Return the generated uuid NOT THREAD SAFE

27.334.2.2 static bool gdcmm::UUIDGenerator::IsValid ( const char \* uid ) [static]

Find out if the string is a valid UUID or not.

The documentation for this class was generated from the following file:


- [gdcmmUUIDGenerator.h](#)

## 27.335 gdcm::Validate Class Reference

[Validate](#) class.

```
#include <gdcmValidate.h>
```

Collaboration diagram for gdcm::Validate:

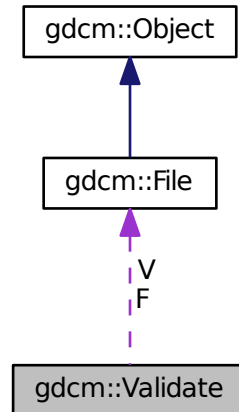

### Public Member Functions

- [Validate](#) ()
- [~Validate](#) ()
- const [File](#) & [GetValidatedFile](#) ()
- void [SetFile](#) ([File](#) const &f)
- void [Validation](#) ()

### Protected Attributes

- const [File](#) \* [F](#)
- [File](#) [V](#)

### 27.335.1 Detailed Description

[Validate](#) class.

### 27.335.2 Constructor & Destructor Documentation

#### 27.335.2.1 gdcm::Validate::Validate ( )

27.335.2.2 `gdcm::Validate::~~Validate ( )`

### 27.335.3 Member Function Documentation

27.335.3.1 `const File& gdcm::Validate::GetValidatedFile ( )` `[inline]`

27.335.3.2 `void gdcm::Validate::SetFile ( File const & f )` `[inline]`

27.335.3.3 `void gdcm::Validate::Validation ( )`

### 27.335.4 Member Data Documentation

27.335.4.1 `const File* gdcm::Validate::F` `[protected]`

27.335.4.2 `File gdcm::Validate::V` `[protected]`

The documentation for this class was generated from the following file:


- [gdcmValidate.h](#)

## 27.336 gdcm::Value Class Reference


Class to represent the value of a Data [Element](#).

```
#include <gdcmValue.h>
```

Inheritance diagram for `gdcm::Value`:

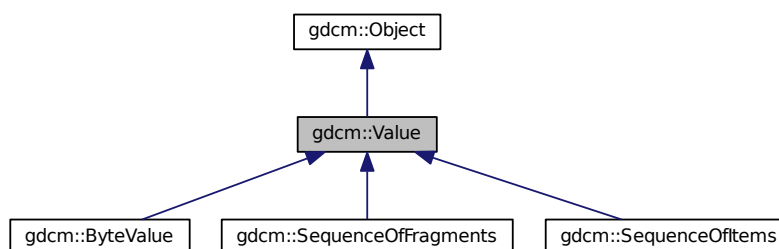

Collaboration diagram for gdcm::Value:

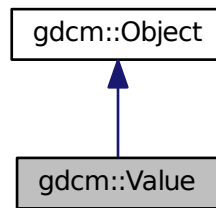

### Public Member Functions

- [Value](#) ()
- [~Value](#) ()
- virtual void [Clear](#) ()=0
- virtual [VL GetLength](#) () const =0
- virtual bool [operator==](#) (const [Value](#) &val) const =0
- virtual void [SetLength](#) ([VL](#) l)=0

### Protected Member Functions

- virtual void [SetLengthOnly](#) ([VL](#) l)

### Friends

- class [DataElement](#)

## 27.336.1 Detailed Description

Class to represent the value of a Data [Element](#).

### Note

VALUE: A component of a [Value](#) Field. A [Value](#) Field may consist of one or more of these components.

## 27.336.2 Constructor & Destructor Documentation

27.336.2.1 `gdcm::Value::Value ( )` `[inline]`

27.336.2.2 `gdcm::Value::~~Value ( )` `[inline]`

## 27.336.3 Member Function Documentation

27.336.3.1 `virtual void gdcm::Value::Clear ( ) [pure virtual]`

Implemented in [gdcm::ByteValue](#), [gdcm::SequenceOfItems](#), and [gdcm::SequenceOfFragments](#).

27.336.3.2 `virtual VL gdcm::Value::GetLength ( ) const [pure virtual]`

Implemented in [gdcm::ByteValue](#), [gdcm::SequenceOfItems](#), and [gdcm::SequenceOfFragments](#).

Referenced by [gdcm::DataSet::InsertDataElement\(\)](#), and [gdcm::DataElement::SetValue\(\)](#).

27.336.3.3 `virtual bool gdcm::Value::operator== ( const Value & val ) const [pure virtual]`

Implemented in [gdcm::SequenceOfFragments](#), [gdcm::SequenceOfItems](#), and [gdcm::ByteValue](#).

27.336.3.4 `virtual void gdcm::Value::SetLength ( VL / ) [pure virtual]`

Implemented in [gdcm::ByteValue](#), [gdcm::SequenceOfItems](#), and [gdcm::SequenceOfFragments](#).

27.336.3.5 `virtual void gdcm::Value::SetLengthOnly ( VL / ) [protected],[virtual]`

Reimplemented in [gdcm::ByteValue](#).

## 27.336.4 Friends And Related Function Documentation

27.336.4.1 `friend class DataElement [friend]`

The documentation for this class was generated from the following file:

- [gdcmValue.h](#)

## 27.337 `gdcm::ValueIO< TDE, TSwap, TType >` Class Template Reference

Class to dispatch template calls.

```
#include <gdcmValueIO.h>
```

### Static Public Member Functions

- static `std::istream & Read (std::istream &is, Value &v, bool readvalues)`
- static `const std::ostream & Write (std::ostream &os, const Value &v)`

### 27.337.1 Detailed Description

```
template<typename TDE, typename TSwap, typename TType = uint8_t>class gdcm::ValueIO< TDE, TSwap, TType >
```

Class to dispatch template calls.

## 27.337.2 Member Function Documentation

27.337.2.1 `template<typename TDE , typename TSwap , typename TType = uint8_t> static std::istream& gdcm::ValueIO< TDE, TSwap, TType >::Read ( std::istream & is, Value & v, bool readvalues ) [static]`

27.337.2.2 `template<typename TDE , typename TSwap , typename TType = uint8_t> static const std::ostream& gdcm::ValueIO< TDE, TSwap, TType >::Write ( std::ostream & os, const Value & v ) [static]`

The documentation for this class was generated from the following file:

- [gdcmValueIO.h](#)

## 27.338 gdcm::Version Class Reference

major/minor and build version

```
#include <gdcmVersion.h>
```

### Public Member Functions

- [Version](#) ()
- [~Version](#) ()
- void [Print](#) (std::ostream &os=std::cout) const

### Static Public Member Functions

- static int [GetBuildVersion](#) ()
- static int [GetMajorVersion](#) ()
- static int [GetMinorVersion](#) ()
- static const char \* [GetVersion](#) ()

### Friends

- std::ostream & [operator<<](#) (std::ostream &\_os, const [Version](#) &v)

## 27.338.1 Detailed Description

major/minor and build version

## 27.338.2 Constructor & Destructor Documentation

27.338.2.1 `gdcm::Version::Version ( ) [inline]`

27.338.2.2 `gdcm::Version::~~Version ( ) [inline]`

## 27.338.3 Member Function Documentation

- 27.338.3.1 `static int gdcmm::Version::GetBuildVersion ( ) [static]`
- 27.338.3.2 `static int gdcmm::Version::GetMajorVersion ( ) [static]`
- 27.338.3.3 `static int gdcmm::Version::GetMinorVersion ( ) [static]`
- 27.338.3.4 `static const char* gdcmm::Version::GetVersion ( ) [static]`
- 27.338.3.5 `void gdcmm::Version::Print ( std::ostream & os = std::cout ) const`

Referenced by `gdcmm::operator<<()`.

## 27.338.4 Friends And Related Function Documentation

- 27.338.4.1 `std::ostream& operator<< ( std::ostream & _os, const Version & v ) [friend]`

The documentation for this class was generated from the following file:

- [gdcmmVersion.h](#)

## 27.339 gdcmm::VL Class Reference

[Value](#) Length.

```
#include <gdcmmVL.h>
```

### Public Types

- `typedef uint32_t Type`

### Public Member Functions

- [VL](#) (uint32\_t vl=0)
- [VL GetLength](#) () const
- `bool IsOdd () const`  
*Return whether or not the [VL](#) is odd or not.*
- `bool IsUndefined () const`
- `operator uint32_t () const`
- [VL & operator++](#) ()
- [VL operator++](#) (int)
- [VL & operator+=](#) ([VL](#) const &vl)  
*+= operator*
- `template<typename TSwap >`  
`std::istream & Read (std::istream &is)`
- `template<typename TSwap >`  
`std::istream & Read16 (std::istream &is)`
- `void SetToUndefined ()`
- `template<typename TSwap >`  
`const std::ostream & Write (std::ostream &os) const`



- `template<typename TSwap >`  
`const std::ostream & Write16 (std::ostream &os) const`

### Static Public Member Functions

- `static uint16_t GetVL16Max ()`
- `static uint32_t GetVL32Max ()`

### Friends

- `std::ostream & operator<< (std::ostream &os, const VL &vl)`

## 27.339.1 Detailed Description

[Value](#) Length.

#### Warning

this is a 4bytes value ! Do not try to use it for 2bytes value length

#### Examples:

[rle2img.cxx](#).

## 27.339.2 Member Typedef Documentation

27.339.2.1 `typedef uint32_t gdcm::VL::Type`

## 27.339.3 Constructor & Destructor Documentation

27.339.3.1 `gdcm::VL::VL ( uint32_t vl = 0 ) \[inline\]`

## 27.339.4 Member Function Documentation

27.339.4.1 `VL gdcm::VL::GetLength ( ) const \[inline\]`

Referenced by `gdcm::FileMetaInformation::GetFullLength()`, and `gdcm::Item::Write()`.

27.339.4.2 `static uint16_t gdcm::VL::GetVL16Max ( ) \[inline\], \[static\]`

27.339.4.3 `static uint32_t gdcm::VL::GetVL32Max ( ) \[inline\], \[static\]`

27.339.4.4 `bool gdcm::VL::IsOdd ( ) const \[inline\]`

Return whether or not the [VL](#) is odd or not.

Referenced by `gdcm::ByteValue::SetLength()`.

27.339.4.5 `bool gdcM::VL::IsUndefined ( ) const [inline]`

Referenced by `gdcM::ByteValue::SetLength()`.

27.339.4.6 `gdcM::VL::operator uint32_t ( ) const [inline]`

27.339.4.7 `VL& gdcM::VL::operator++ ( ) [inline]`

27.339.4.8 `VL gdcM::VL::operator++ ( int ) [inline]`

27.339.4.9 `VL& gdcM::VL::operator+= ( VL const & v/ ) [inline]`

`+=` operator

27.339.4.10 `template<typename TSwap > std::istream& gdcM::VL::Read ( std::istream & is ) [inline]`

27.339.4.11 `template<typename TSwap > std::istream& gdcM::VL::Read16 ( std::istream & is ) [inline]`

27.339.4.12 `void gdcM::VL::SetToUndefined ( ) [inline]`

27.339.4.13 `template<typename TSwap > const std::ostream& gdcM::VL::Write ( std::ostream & os ) const [inline]`

Referenced by `gdcM::Fragment::Write()`, `gdcM::SequenceOfItems::Write()`, `gdcM::Item::Write()`, and `gdcM::SequenceOfFragments::Write()`.

27.339.4.14 `template<typename TSwap > const std::ostream& gdcM::VL::Write16 ( std::ostream & os ) const [inline]`

## 27.339.5 Friends And Related Function Documentation

27.339.5.1 `std::ostream& operator<< ( std::ostream & os, const VL & v/ ) [friend]`

The documentation for this class was generated from the following file:

- [gdcMVL.h](#)

## 27.340 gdcM::VM Class Reference

**Value** Multiplicity Looking at the DICOMV3 dict only there is very few cases: 1 2 3 4 5 6 8 16 24 1-2 1-3 1-8 1-32 1-99 1-n 2-2n 2-n 3-3n 3-n.

```
#include <gdcMVM.h>
```

## Public Types

- enum [VMType](#) {
  - [VM0](#) = 0,
  - [VM1](#) = 1,
  - [VM2](#) = 2,
  - [VM3](#) = 4,
  - [VM4](#) = 8,
  - [VM5](#) = 16,
  - [VM6](#) = 32,
  - [VM8](#) = 64,
  - [VM9](#) = 128,
  - [VM10](#) = 256,
  - [VM12](#) = 512,
  - [VM16](#) = 1024,
  - [VM18](#) = 2048,
  - [VM24](#) = 4096,
  - [VM28](#) = 8192,
  - [VM32](#) = 16384,
  - [VM35](#) = 32768,
  - [VM99](#) = 65536,
  - [VM256](#) = 131072,
  - [VM1\\_2](#) = VM1 | VM2,
  - [VM1\\_3](#) = VM1 | VM2 | VM3,
  - [VM1\\_4](#) = VM1 | VM2 | VM3 | VM4,
  - [VM1\\_5](#) = VM1 | VM2 | VM3 | VM4 | VM5,
  - [VM1\\_8](#) = VM1 | VM2 | VM3 | VM4 | VM5 | VM6 | VM8,
  - [VM1\\_32](#) = VM1 | VM2 | VM3 | VM4 | VM5 | VM6 | VM8 | VM9 | VM16 | VM24 | VM32,
  - [VM1\\_99](#) = VM1 | VM2 | VM3 | VM4 | VM5 | VM6 | VM8 | VM9 | VM16 | VM24 | VM32 | VM99,
  - [VM1\\_n](#) = VM1 | VM2 | VM3 | VM4 | VM5 | VM6 | VM8 | VM9 | VM16 | VM24 | VM32 | VM99 | VM256,
  - [VM2\\_2n](#) = VM2 | VM4 | VM6 | VM8 | VM16 | VM24 | VM32 | VM256,
  - [VM2\\_n](#) = VM2 | VM3 | VM4 | VM5 | VM6 | VM8 | VM9 | VM16 | VM24 | VM32 | VM99 | VM256,
  - [VM3\\_4](#) = VM3 | VM4,
  - [VM3\\_3n](#) = VM3 | VM6 | VM9 | VM24 | VM99 | VM256,
  - [VM3\\_n](#) = VM3 | VM4 | VM5 | VM6 | VM8 | VM9 | VM16 | VM24 | VM32 | VM99 | VM256,
  - [VM4\\_4n](#) = VM4 | VM16 | VM24 | VM32 | VM256,
  - [VM6\\_6n](#) = VM6 | VM12 | VM18 | VM24,
  - [VM7\\_7n](#),
  - [VM30\\_30n](#),
  - [VM47\\_47n](#),
  - [VM\\_END](#) = VM1\_n + 1 }

## Public Member Functions

- [VM](#) ([VMType](#) type=[VM0](#))
- bool [Compatible](#) ([VM](#) const &vm) const
- unsigned int [GetLength](#) () const
- operator [VMType](#) () const

## Static Public Member Functions

- static unsigned int [GetNumberOfElementsFromArray](#) (const char \*array, unsigned int length)

- static const char \* [GetVMString](#) ([VMType](#) vm)
- static [VMType](#) [GetVMType](#) (const char \*vm)
- static [VMType](#) [GetVMTypeFromLength](#) (unsigned int length, unsigned int size)
- static bool [IsValid](#) (int vm1, [VMType](#) vm2)

### Static Protected Member Functions

- static unsigned int [GetIndex](#) ([VMType](#) vm)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [VM](#) &vm)

### 27.340.1 Detailed Description

[Value](#) Multiplicity Looking at the DICOMV3 dict only there is very few cases: 1 2 3 4 5 6 8 16 24 1-2 1-3 1-8 1-32 1-99 1-n 2-2n 2-n 3-3n 3-n.

Some private dict define some more: 4-4n 1-4 1-5 256 9 3-4

even more:

7-7n 10 18 12 35 47\_47n 30\_30n 28

6-6n

### 27.340.2 Member Enumeration Documentation

#### 27.340.2.1 enum `gdcm::VM::VMType`

Enumerator

***VM0***  
***VM1***  
***VM2***  
***VM3***  
***VM4***  
***VM5***  
***VM6***  
***VM8***  
***VM9***  
***VM10***  
***VM12***  
***VM16***  
***VM18***  
***VM24***  
***VM28***  
***VM32***

**VM35**  
**VM99**  
**VM256**  
**VM1\_2**  
**VM1\_3**  
**VM1\_4**  
**VM1\_5**  
**VM1\_8**  
**VM1\_32**  
**VM1\_99**  
**VM1\_n**  
**VM2\_2n**  
**VM2\_n**  
**VM3\_4**  
**VM3\_3n**  
**VM3\_n**  
**VM4\_4n**  
**VM6\_6n**  
**VM7\_7n**  
**VM30\_30n**  
**VM47\_47n**  
**VM\_END**

### 27.340.3 Constructor & Destructor Documentation

27.340.3.1 `gdcmm::VM ( VMType type = VM0 ) [inline]`

### 27.340.4 Member Function Documentation

27.340.4.1 `bool gdcmm::VM::Compatible ( VM const & vm ) const`

WARNING: Implementation deficiency The Compatible function is poorly implemented, the reference vm should be coming from the dictionary, while the passed in value is the value guess from the file.

27.340.4.2 `static unsigned int gdcmm::VM::GetIndex ( VMType vm ) [static], [protected]`

27.340.4.3 `unsigned int gdcmm::VM::GetLength ( ) const`

27.340.4.4 `static unsigned int gdcmm::VM::GetNumberOfElementsFromArray ( const char * array, unsigned int length ) [static]`

27.340.4.5 `static const char* gdcmm::VM::GetVMString ( VMType vm ) [static]`

Return the string as written in the official DICOM dict from a custom enum type

Referenced by `gdcmm::operator<<()`.

27.340.4.6 static VMType gdcM::VM::GetVMType ( const char \* *vm* ) [static]

27.340.4.7 static VMType gdcM::VM::GetVMTypeFromLength ( unsigned int *length*, unsigned int *size* ) [static]

27.340.4.8 static bool gdcM::VM::IsValid ( int *vm1*, VMType *vm2* ) [static]

Check if *vm1* is valid compare to *vm2*, i.e *vm1* is element of *vm2* *vm1* is typically deduce from counting in a ValueField

27.340.4.9 gdcM::VM::operator VMType ( ) const [inline]

## 27.340.5 Friends And Related Function Documentation

27.340.5.1 std::ostream& operator<< ( std::ostream & *os*, const VM & *vm* ) [friend]

The documentation for this class was generated from the following file:

- [gdcMVM.h](#)

## 27.341 gdcM::VMToLength< T > Struct Template Reference

```
#include <gdcMVM.h>
```

The documentation for this struct was generated from the following file:

- [gdcMVM.h](#)

## 27.342 gdcM::VR Class Reference

[VR](#) class This is adapted from DICOM standard The biggest difference is the INVALID [VR](#) and the composite one that differ from standard (more like an addition) This allow us to represent all the possible case express in the DICOMV3 dict.

```
#include <gdcMVR.h>
```

## Public Types

- enum [VRType](#) {
  - [INVALID](#) = 0,
  - [AE](#) = 1,
  - [AS](#) = 2,
  - [AT](#) = 4,
  - [CS](#) = 8,
  - [DA](#) = 16,
  - [DS](#) = 32,
  - [DT](#) = 64,
  - [FD](#) = 128,
  - [FL](#) = 256,
  - [IS](#) = 512,
  - [LO](#) = 1024,
  - [LT](#) = 2048,
  - [OB](#) = 4096,
  - [OD](#) = 134217728,
  - [OF](#) = 8192,
  - [OW](#) = 16384,
  - [PN](#) = 32768,
  - [SH](#) = 65536,
  - [SL](#) = 131072,
  - [SQ](#) = 262144,
  - [SS](#) = 524288,
  - [ST](#) = 1048576,
  - [TM](#) = 2097152,
  - [UI](#) = 4194304,
  - [UL](#) = 8388608,
  - [UN](#) = 16777216,
  - [US](#) = 33554432,
  - [UT](#) = 67108864,
  - [OB\\_OW](#) = OB | OW,
  - [US\\_SS](#) = US | SS,
  - [US\\_SS\\_OW](#) = US | SS | OW,
  - [VL16](#) = AE | AS | AT | CS | DA | DS | DT | FD | FL | IS | LO | LT | PN | SH | SL | SS | ST | TM | UI | UL | US,
  - [VL32](#) = OB | OW | OD | OF | SQ | UN | UT,
  - [VRASCII](#) = AE | AS | CS | DA | DS | DT | IS | LO | LT | PN | SH | ST | TM | UI | UT,
  - [VRBINARY](#) = AT | FL | FD | OB | OD | OF | OW | SL | SQ | SS | UL | UN | US,
  - [VR\\_VM1](#) = AS | LT | ST | UT | SQ | OF | OD | OW | OB | UN,
  - [VRALL](#) = VRASCII | VRBINARY,
  - [VR\\_END](#) = UT+1 }

## Public Member Functions

- [VR](#) ([VRType](#) vr=[INVALID](#))
- bool [Compatible](#) ([VR](#) const &vr) const
- int [GetLength](#) () const
- unsigned int [GetSize](#) () const
- unsigned int [GetSizeof](#) () const
- bool [IsDual](#) () const
- bool [IsVRFile](#) () const
- [operator VRType](#) () const

- std::istream & [Read](#) (std::istream &is)
- const std::ostream & [Write](#) (std::ostream &os) const

### Static Public Member Functions

- static bool [CanDisplay](#) (VRType vr)
- static uint32\_t [GetLength](#) (VRType vr)
- static const char \* [GetVRString](#) (VRType vr)
- static const char \* [GetVRStringFromFile](#) (VRType vr)
- static VRType [GetVRType](#) (const char \*vr)
- static VRType [GetVRTypeFromFile](#) (const char \*vr)
- static bool [IsASCII](#) (VRType vr)
- static bool [IsASCII2](#) (VRType vr)
- static bool [IsBinary](#) (VRType vr)
- static bool [IsBinary2](#) (VRType vr)
- static bool [IsSwap](#) (const char \*vr)
- static bool [IsValid](#) (const char \*vr)
- static bool [IsValid](#) (const char \*vr1, VRType vr2)

### Friends

- std::ostream & [operator<<](#) (std::ostream &os, const [VR](#) &vr)

## 27.342.1 Detailed Description

[VR](#) class This is adapted from DICOM standard The biggest difference is the INVALID [VR](#) and the composite one that differ from standard (more like an addition) This allow us to represent all the possible case express in the DICOMV3 dict.

#### Note

VALUE REPRESENTATION ([VR](#)) Specifies the data type and format of the Value(s) contained in the [Value](#) Field of a Data [Element](#). VALUE REPRESENTATION FIELD: The field where the [Value](#) Representation of a Data [Element](#) is stored in the encoding of a Data [Element](#) structure with explicit [VR](#).

#### Examples:

[GenAllVR.cxx](#), and [GenFakeIdentifyFile.cxx](#).

## 27.342.2 Member Enumeration Documentation

### 27.342.2.1 enum gdcm::VR::VRType

#### Enumerator

**INVALID**  
**AE**  
**AS**  
**AT**  
**CS**



*DA*  
*DS*  
*DT*  
*FD*  
*FL*  
*IS*  
*LO*  
*LT*  
*OB*  
*OD*  
*OF*  
*OW*  
*PN*  
*SH*  
*SL*  
*SQ*  
*SS*  
*ST*  
*TM*  
*UI*  
*UL*  
*UN*  
*US*  
*UT*  
*OB\_OW*  
*US\_SS*  
*US\_SS\_OW*  
*VL16*  
*VL32*  
*VRASCII*  
*VRBINARY*  
*VR\_VM1*  
*VRALL*  
*VR\_END*

Examples:

[NewSequence.cs](#).

### 27.342.3 Constructor & Destructor Documentation

27.342.3.1 `gdcm::VR::VR ( VRType vr = INVALID ) [inline]`

### 27.342.4 Member Function Documentation

27.342.4.1 `static bool gdcm::VR::CanDisplay ( VRType vr ) [static]`

27.342.4.2 `bool gdcm::VR::Compatible ( VR const & vr ) const`

27.342.4.3 `int gdcm::VR::GetLength ( ) const [inline]`

27.342.4.4 `static uint32_t gdcm::VR::GetLength ( VRType vr ) [inline],[static]`

27.342.4.5 `unsigned int gdcm::VR::GetSize ( ) const [inline]`

References AE, US\_SS, and VRTypeTemplateCase.

27.342.4.6 `unsigned int gdcm::VR::GetSizeof ( ) const`

27.342.4.7 `static const char* gdcm::VR::GetVRString ( VRType vr ) [static]`

Referenced by `gdcm::operator<<()`.

27.342.4.8 `static const char* gdcm::VR::GetVRStringFromFile ( VRType vr ) [static]`

27.342.4.9 `static VRType gdcm::VR::GetVRType ( const char * vr ) [static]`

27.342.4.10 `static VRType gdcm::VR::GetVRTypeFromFile ( const char * vr ) [static]`

27.342.4.11 `static bool gdcm::VR::IsASCII ( VRType vr ) [static]`

27.342.4.12 `static bool gdcm::VR::IsASCII2 ( VRType vr ) [static]`

27.342.4.13 `static bool gdcm::VR::IsBinary ( VRType vr ) [static]`

27.342.4.14 `static bool gdcm::VR::IsBinary2 ( VRType vr ) [static]`

27.342.4.15 `bool gdcm::VR::IsDual ( ) const`

27.342.4.16 `static bool gdcm::VR::IsSwap ( const char * vr ) [static]`

27.342.4.17 `static bool gdcm::VR::IsValid ( const char * vr ) [static]`

27.342.4.18 `static bool gdcm::VR::IsValid ( const char * vr1, VRType vr2 ) [static]`

27.342.4.19 `bool gdcm::VR::IsVRFile ( ) const`

Referenced by `gdcm::DataElement::SetVR()`.

27.342.4.20 `gdcm::VR::operator VRType ( ) const` `[inline]`

27.342.4.21 `std::istream& gdcm::VR::Read ( std::istream & is )` `[inline]`

References `gdcmDebugMacro`, `INVALID`, and `VR_END`.

27.342.4.22 `const std::ostream& gdcm::VR::Write ( std::ostream & os ) const` `[inline]`

References `gdcmAssertAlwaysMacro`, and `INVALID`.

## 27.342.5 Friends And Related Function Documentation

27.342.5.1 `std::ostream& operator<< ( std::ostream & os, const VR & vr )` `[friend]`

The documentation for this class was generated from the following file:


- [gdcmVR.h](#)

## 27.343 gdcm::VR16ExplicitDataElement Class Reference


Class to read/write a [DataElement](#) as Explicit Data [Element](#).

```
#include <gdcmVR16ExplicitDataElement.h>
```

Inheritance diagram for `gdcm::VR16ExplicitDataElement`:

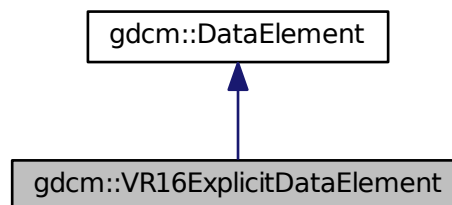

Collaboration diagram for `gdcm::VR16ExplicitDataElement`:

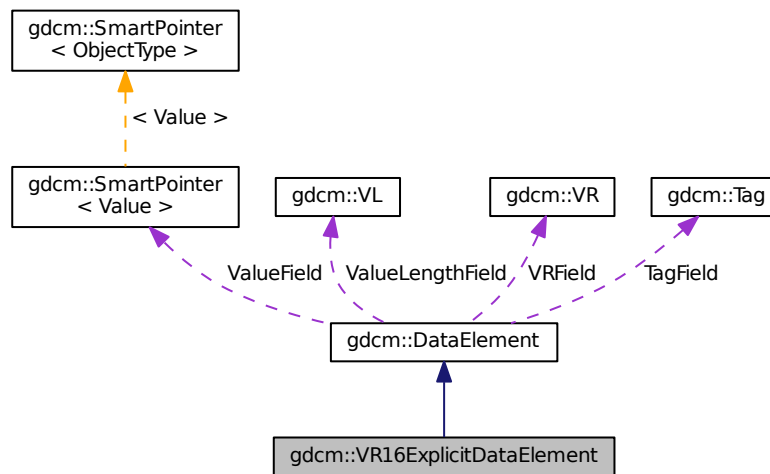

## Public Member Functions

- [VL GetLength](#) () const
- template<typename TSwap >  
std::istream & [Read](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadPreValue](#) (std::istream &is)
- template<typename TSwap >  
std::istream & [ReadValue](#) (std::istream &is, bool readvalues=true)
- template<typename TSwap >  
std::istream & [ReadWithLength](#) (std::istream &is, [VL](#) &length)

## Additional Inherited Members

### 27.343.1 Detailed Description

Class to read/write a [DataElement](#) as Explicit Data [Element](#).

#### Note

This class support 16 bits when finding an unkown [VR](#): For instance: Siemens\_CT\_Sensation64\_has\_VR\_RT.dcm

### 27.343.2 Member Function Documentation

27.343.2.1 [VL](#) `gdcm::VR16ExplicitDataElement::GetLength ( )` const

27.343.2.2 template<typename TSwap > std::istream& `gdcm::VR16ExplicitDataElement::Read ( std::istream & is )`

27.343.2.3 `template<typename TSwap > std::istream& gdcm::VR16ExplicitDataElement::ReadPreValue ( std::istream & is )`

27.343.2.4 `template<typename TSwap > std::istream& gdcm::VR16ExplicitDataElement::ReadValue ( std::istream & is, bool readvalues = true )`

27.343.2.5 `template<typename TSwap > std::istream& gdcm::VR16ExplicitDataElement::ReadWithLength ( std::istream & is, VL & length )`

The documentation for this class was generated from the following file:

- [gdcmVR16ExplicitDataElement.h](#)

## 27.344 gdcm::VRToEncoding< T > Struct Template Reference

```
#include <gdcmVR.h>
```


The documentation for this struct was generated from the following file:

- [gdcmVR.h](#)

## 27.345 gdcm::VRToType< T > Struct Template Reference

```
#include <gdcmVR.h>
```

Inheritance diagram for gdcm::VRToType< T >:

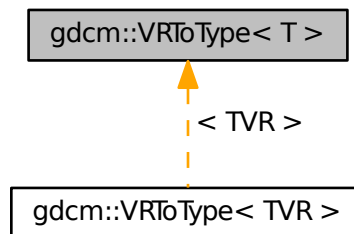

### 27.345.1 Detailed Description

```
template<int T>struct gdcm::VRToType< T >
```

Examples:

[DumpGEMSMovieGroup.cxx](#).

The documentation for this struct was generated from the following file:

- [gdcmVR.h](#)

## 27.346 `gdcm::VRVLSIZE< T >` Class Template Reference

```
#include <gdcmAttribute.h>
```

The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)

## 27.347 `gdcm::VRVLSIZE< 0 >` Class Template Reference

```
#include <gdcmAttribute.h>
```

### Static Public Member Functions

- static `uint16_t` [Read](#) (`std::istream &_is`)
- static void [Write](#) (`std::ostream &os`)

### 27.347.1 Member Function Documentation

27.347.1.1 static `uint16_t` `gdcm::VRVLSIZE< 0 >::Read` (`std::istream &_is`) [`inline`], [`static`]

27.347.1.2 static void `gdcm::VRVLSIZE< 0 >::Write` (`std::ostream &os`) [`inline`], [`static`]

The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)

## 27.348 `gdcm::VRVLSIZE< 1 >` Class Template Reference

```
#include <gdcmAttribute.h>
```

### Static Public Member Functions

- static `uint32_t` [Read](#) (`std::istream &_is`)
- static void [Write](#) (`std::ostream &os`)

### 27.348.1 Member Function Documentation

27.348.1.1 static `uint32_t` `gdcm::VRVLSIZE< 1 >::Read` (`std::istream &_is`) [`inline`], [`static`]

27.348.1.2 static void `gdcm::VRVLSIZE< 1 >::Write` (`std::ostream &os`) [`inline`], [`static`]


The documentation for this class was generated from the following file:

- [gdcmAttribute.h](#)


## 27.349 vtkGDCMImageReader Class Reference

```
#include <vtkGDCMImageReader.h>
```

Inheritance diagram for vtkGDCMImageReader:

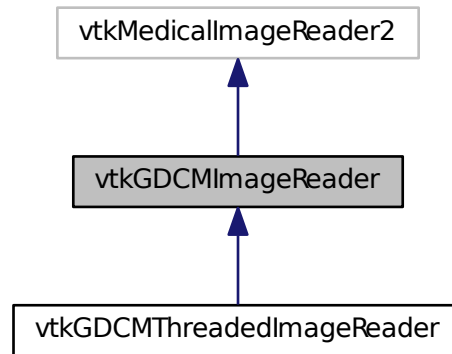

Collaboration diagram for vtkGDCMImageReader:

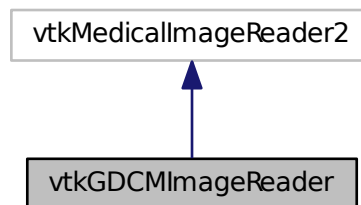

### Public Member Functions

- virtual int [CanReadFile](#) (const char \*fname)
- virtual const char \* [GetDescriptiveName](#) ()
- virtual const char \* [GetFileExtensions](#) ()
- vtkImageData \* [GetIconImage](#) ()
- vtkImageData \* [GetOverlay](#) (int i)
- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetCurve](#) (vtkPolyData \*pd)
- virtual void [SetFileNames](#) (vtkStringArray \*)

- virtual void [SetMedicalImageProperties](#) (vtkMedicalImageProperties \*pd)
- [vtkBooleanMacro](#) (LoadOverlays, int)
- [vtkBooleanMacro](#) (LoadIconImage, int)
- [vtkBooleanMacro](#) (LossyFlag, int)
- [vtkBooleanMacro](#) (ApplyLookupTable, int)
- int [vtkBooleanMacro](#) (ApplyYBRToRGB, int)
- [vtkGetMacro](#) (LoadOverlays, int)
- [vtkGetMacro](#) (LoadIconImage, int)
- [vtkGetMacro](#) (LossyFlag, int)
- [vtkGetMacro](#) (NumberOfOverlays, int)
- [vtkGetMacro](#) (NumberOfIconImages, int)
- [vtkGetMacro](#) (ApplyLookupTable, int)
- [vtkGetMacro](#) (ApplyYBRToRGB, int) [vtkSetMacro](#)(ApplyYBRToRGB
- [vtkGetMacro](#) (ImageFormat, int)
- [vtkGetMacro](#) (PlanarConfiguration, int)
- [vtkGetMacro](#) (Shift, double)
- [vtkGetMacro](#) (Scale, double)
- [vtkGetObjectMacro](#) (DirectionCosines, vtkMatrix4x4)
- [vtkGetObjectMacro](#) (MedicalImageProperties, vtkMedicalImageProperties)
- [vtkGetObjectMacro](#) (FileNames, vtkStringArray)
- [vtkGetObjectMacro](#) (Curve, vtkPolyData)
- [vtkGetVector3Macro](#) (ImagePositionPatient, double)
- [vtkGetVector6Macro](#) (ImageOrientationPatient, double)
- [vtkSetMacro](#) (LoadOverlays, int)
- [vtkSetMacro](#) (LoadIconImage, int)
- [vtkSetMacro](#) (LossyFlag, int)
- [vtkSetMacro](#) (ApplyLookupTable, int)
- [vtkTypeRevisionMacro](#) (vtkGDCMImageReader, vtkMedicalImageReader2)

### Static Public Member Functions

- static [vtkGDCMImageReader](#) \* [New](#) ()

### Protected Member Functions

- [vtkGDCMImageReader](#) ()
- [~vtkGDCMImageReader](#) ()
- void [ExecuteData](#) (vtkDataObject \*out)
- void [ExecuteInformation](#) ()
- void [FillMedicalImageInformation](#) (const [gdcm::ImageReader](#) &reader)
- int [LoadSingleFile](#) (const char \*filename, char \*pointer, unsigned long &outlen)
- int [RequestDataCompat](#) ()
- int [RequestInformationCompat](#) ()
- void [SetFilePattern](#) (const char \*)
- void [SetFilePrefix](#) (const char \*)
- [vtkGetStringMacro](#) (FilePrefix)
- [vtkGetStringMacro](#) (FilePattern)
- [vtkSetVector6Macro](#) (ImageOrientationPatient, double)



## Protected Attributes

- int [ApplyInverseVideo](#)
- int [ApplyLookupTable](#)
- int [ApplyPlanarConfiguration](#)
- int [ApplyShiftScale](#)
- int [ApplyYBRToRGB](#)
- vtkPolyData \* [Curve](#)
- vtkMatrix4x4 \* [DirectionCosines](#)
- vtkStringArray \* [FileNames](#)
- int [ForceRescale](#)
- int [IconDataScalarType](#)
- int [IconImageDataExtent](#) [6]
- int [IconNumberOfScalarComponents](#)
- int [ImageFormat](#)
- double [ImageOrientationPatient](#) [6]
- double [ImagePositionPatient](#) [3]
- int [LoadIconImage](#)
- int [LoadOverlays](#)
- int [LossyFlag](#)
- vtkMedicalImageProperties \* [MedicalImageProperties](#)
- int [NumberOfIconImages](#)
- int [NumberOfOverlays](#)
- int [PlanarConfiguration](#)
- double [Scale](#)
- double [Shift](#)

### 27.349.1 Detailed Description

Examples:

[AWTMedical3.java](#), [Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [gdcmmorthoplanes.cxx](#), [gdcmmreslice.cxx](#), [gdcmmtexture.cxx](#), [gdcmmvolume.cxx](#), [HelloActiviz.cs](#), [HelloActiviz2.cs](#), [HelloActiviz3.cs](#), [HelloActiviz4.cs](#), [HelloActiviz5.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld.java](#), [MagnifyFile.cxx](#), [MetaImageMD5Activiz.cs](#), [MIPViewer.java](#), [MPRViewer.java](#), [MPRViewer2.java](#), [offscreenimage.cxx](#), [ReadSeriesIntoVTK.java](#), [RefCounting.cs](#), and [reslicesphere.cxx](#).

### 27.349.2 Constructor & Destructor Documentation

27.349.2.1 `vtkGDCMImageReader::vtkGDCMImageReader ( )` [protected]

Examples:

[HelloActiviz2.cs](#).

27.349.2.2 `vtkGDCMImageReader::~~vtkGDCMImageReader ( )` [protected]

### 27.349.3 Member Function Documentation

27.349.3.1 `virtual int vtkGDCMImageReader::CanReadFile ( const char * fname )` [virtual]

Examples:

[MetalImageMD5Activiz.cs](#).

27.349.3.2 `void vtkGDCMImageReader::ExecuteData ( vtkDataObject * out )` [protected]

27.349.3.3 `void vtkGDCMImageReader::ExecuteInformation ( )` [protected]

27.349.3.4 `void vtkGDCMImageReader::FillMedicalImageInformation ( const gdcm::ImageReader & reader )`  
[protected]

27.349.3.5 `virtual const char* vtkGDCMImageReader::GetDescriptiveName ( )` [inline],[virtual]

27.349.3.6 `virtual const char* vtkGDCMImageReader::GetFileExtensions ( )` [inline],[virtual]

27.349.3.7 `vtkImageData* vtkGDCMImageReader::GetIconImage ( )`

27.349.3.8 `vtkImageData* vtkGDCMImageReader::GetOverlay ( int i )`

27.349.3.9 `int vtkGDCMImageReader::LoadSingleFile ( const char * filename, char * pointer, unsigned long & outlen )`  
[protected]

27.349.3.10 `static vtkGDCMImageReader* vtkGDCMImageReader::New ( )` [static]

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [gdcmmorphoplanes.cxx](#), [gdcmreslice.cxx](#), [gdcmttexture.cxx](#), [gdcmvolume.cxx](#), [HelloActiviz.cs](#), [HelloActiviz3.cs](#), [HelloActiviz4.cs](#), [HelloActiviz5.cs](#), [HelloVTKWorld.cs](#), [MagnifyFile.cxx](#), [MetalImageMD5Activiz.cs](#), [offscreenimage.cxx](#), [RefCounting.cs](#), and [reslicesphere.cxx](#).

27.349.3.11 `virtual void vtkGDCMImageReader::PrintSelf ( ostream & os, vtkIndent indent )` [virtual]

Reimplemented in [vtkGDCMThreadedImageReader](#).

27.349.3.12 `int vtkGDCMImageReader::RequestDataCompat ( )` [protected]

27.349.3.13 `int vtkGDCMImageReader::RequestInformationCompat ( )` [protected]

27.349.3.14 `virtual void vtkGDCMImageReader::SetCurve ( vtkPolyData * pd )` [virtual]

27.349.3.15 `virtual void vtkGDCMImageReader::SetFileNames ( vtkStringArray * )` [virtual]

Examples:

[gdcmortoplanes.cxx](#), [HelloActiviz3.cs](#), [HelloActiviz4.cs](#), [HelloActiviz5.cs](#), [MIPViewer.java](#), [MPRViewer.java](#), [MPRViewer2.java](#), and [ReadSeriesIntoVTK.java](#).

27.349.3.16 `void vtkGDCMImageReader::SetFilePattern ( const char * )` [inline],[protected]

27.349.3.17 `void vtkGDCMImageReader::SetFilePrefix ( const char * )` [inline],[protected]

27.349.3.18 `virtual void vtkGDCMImageReader::SetMedicalImageProperties ( vtkMedicalImageProperties * pd )` [virtual]

27.349.3.19 `vtkGDCMImageReader::vtkBooleanMacro ( LoadOverlays , int )`

27.349.3.20 `vtkGDCMImageReader::vtkBooleanMacro ( LoadIconImage , int )`

27.349.3.21 `vtkGDCMImageReader::vtkBooleanMacro ( LossyFlag , int )`

27.349.3.22 `vtkGDCMImageReader::vtkBooleanMacro ( ApplyLookupTable , int )`

27.349.3.23 `int vtkGDCMImageReader::vtkBooleanMacro ( ApplyYBRToRGB , int )`

27.349.3.24 `vtkGDCMImageReader::vtkGetMacro ( LoadOverlays , int )`

27.349.3.25 `vtkGDCMImageReader::vtkGetMacro ( LoadIconImage , int )`

27.349.3.26 `vtkGDCMImageReader::vtkGetMacro ( LossyFlag , int )`

27.349.3.27 `vtkGDCMImageReader::vtkGetMacro ( NumberOfOverlays , int )`

27.349.3.28 `vtkGDCMImageReader::vtkGetMacro ( NumberOfIconImages , int )`

27.349.3.29 `vtkGDCMImageReader::vtkGetMacro ( ApplyLookupTable , int )`

27.349.3.30 `vtkGDCMImageReader::vtkGetMacro ( ApplyYBRToRGB , int )`

27.349.3.31 `vtkGDCMImageReader::vtkGetMacro ( ImageFormat , int )`

27.349.3.32 `vtkGDCMImageReader::vtkGetMacro ( PlanarConfiguration , int )`

27.349.3.33 `vtkGDCMImageReader::vtkGetMacro ( Shift , double )`

27.349.3.34 `vtkGDCMImageReader::vtkGetMacro ( Scale , double )`

27.349.3.35 `vtkGDCMImageReader::vtkGetObjectMacro ( DirectionCosines , vtkMatrix4x4 )`

27.349.3.36 `vtkGDCMImageReader::vtkGetObjectMacro ( MedicalImageProperties , vtkMedicalImageProperties )`

27.349.3.37 `vtkGDCMImageReader::vtkGetObjectMacro ( FileNames , vtkStringArray )`

- 27.349.3.38 `vtkGDCMImageReader::vtkGetObjectMacro ( Curve , vtkPolyData )`
- 27.349.3.39 `vtkGDCMImageReader::vtkGetStringMacro ( FilePrefix )` [protected]
- 27.349.3.40 `vtkGDCMImageReader::vtkGetStringMacro ( FilePattern )` [protected]
- 27.349.3.41 `vtkGDCMImageReader::vtkGetVector3Macro ( ImagePositionPatient , double )`
- 27.349.3.42 `vtkGDCMImageReader::vtkGetVector6Macro ( ImageOrientationPatient , double )`
- 27.349.3.43 `vtkGDCMImageReader::vtkSetMacro ( LoadOverlays , int )`
- 27.349.3.44 `vtkGDCMImageReader::vtkSetMacro ( LoadIconImage , int )`
- 27.349.3.45 `vtkGDCMImageReader::vtkSetMacro ( LossyFlag , int )`
- 27.349.3.46 `vtkGDCMImageReader::vtkSetMacro ( ApplyLookupTable , int )`
- 27.349.3.47 `vtkGDCMImageReader::vtkSetVector6Macro ( ImageOrientationPatient , double )` [protected]
- 27.349.3.48 `vtkGDCMImageReader::vtkTypeRevisionMacro ( vtkGDCMImageReader , vtkMedicalImageReader2 )`

#### 27.349.4 Member Data Documentation

- 27.349.4.1 `int vtkGDCMImageReader::ApplyInverseVideo` [protected]
- 27.349.4.2 `int vtkGDCMImageReader::ApplyLookupTable` [protected]
- 27.349.4.3 `int vtkGDCMImageReader::ApplyPlanarConfiguration` [protected]
- 27.349.4.4 `int vtkGDCMImageReader::ApplyShiftScale` [protected]
- 27.349.4.5 `int vtkGDCMImageReader::ApplyYBRToRGB` [protected]
- 27.349.4.6 `vtkPolyData* vtkGDCMImageReader::Curve` [protected]
- 27.349.4.7 `vtkMatrix4x4* vtkGDCMImageReader::DirectionCosines` [protected]
- 27.349.4.8 `vtkStringArray* vtkGDCMImageReader::FileNames` [protected]
- 27.349.4.9 `int vtkGDCMImageReader::ForceRescale` [protected]
- 27.349.4.10 `int vtkGDCMImageReader::IconDataScalarType` [protected]
- 27.349.4.11 `int vtkGDCMImageReader::IconImageDataExtent[6]` [protected]
- 27.349.4.12 `int vtkGDCMImageReader::IconNumberOfScalarComponents` [protected]
- 27.349.4.13 `int vtkGDCMImageReader::ImageFormat` [protected]
- 27.349.4.14 `double vtkGDCMImageReader::ImageOrientationPatient[6]` [protected]

- 27.349.4.15 `double vtkGDCMImageReader::ImagePositionPatient[3]` [protected]
- 27.349.4.16 `int vtkGDCMImageReader::LoadIconImage` [protected]
- 27.349.4.17 `int vtkGDCMImageReader::LoadOverlays` [protected]
- 27.349.4.18 `int vtkGDCMImageReader::LossyFlag` [protected]
- 27.349.4.19 `vtkMedicalImageProperties* vtkGDCMImageReader::MedicalImageProperties` [protected]
- 27.349.4.20 `int vtkGDCMImageReader::NumberOfIconImages` [protected]
- 27.349.4.21 `int vtkGDCMImageReader::NumberOfOverlays` [protected]
- 27.349.4.22 `int vtkGDCMImageReader::PlanarConfiguration` [protected]
- 27.349.4.23 `double vtkGDCMImageReader::Scale` [protected]
- 27.349.4.24 `double vtkGDCMImageReader::Shift` [protected]


The documentation for this class was generated from the following file:

- [vtkGDCMImageReader.h](#)


## 27.350 vtkGDCMImageReader2 Class Reference

```
#include <vtkGDCMImageReader2.h>
```

Inheritance diagram for vtkGDCMImageReader2:

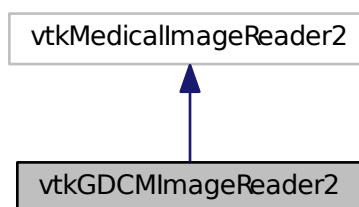

Collaboration diagram for vtkGDCMImageReader2:

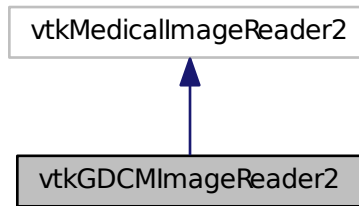

## Public Member Functions

- virtual int [CanReadFile](#) (const char \*fname)
- virtual const char \* [GetDescriptiveName](#) ()
- virtual const char \* [GetFileExtensions](#) ()
- vtkImageData \* [GetIconImage](#) ()
- vtkAlgorithmOutput \* [GetIconImagePort](#) ()
- vtkImageData \* [GetOverlay](#) (int i)
- vtkAlgorithmOutput \* [GetOverlayPort](#) (int index)
- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetCurve](#) (vtkPolyData \*pd)
- virtual void [SetMedicalImageProperties](#) (vtkMedicalImageProperties \*pd)
- [vtkBooleanMacro](#) (LoadOverlays, int)
- [vtkBooleanMacro](#) (LoadIconImage, int)
- [vtkBooleanMacro](#) (LossyFlag, int)
- [vtkBooleanMacro](#) (ApplyLookupTable, int)
- int [vtkBooleanMacro](#) (ApplyYBRToRGB, int)
- [vtkGetMacro](#) (LoadOverlays, int)
- [vtkGetMacro](#) (LoadIconImage, int)
- [vtkGetMacro](#) (LossyFlag, int)
- [vtkGetMacro](#) (NumberOfOverlays, int)
- [vtkGetMacro](#) (NumberOfIconImages, int)
- [vtkGetMacro](#) (ApplyLookupTable, int)
- [vtkGetMacro](#) (ApplyYBRToRGB, int) [vtkSetMacro](#) (ApplyYBRToRGB
- [vtkGetMacro](#) (ImageFormat, int)
- [vtkGetMacro](#) (PlanarConfiguration, int)
- [vtkGetMacro](#) (Shift, double)
- [vtkGetMacro](#) (Scale, double)
- [vtkGetObjectMacro](#) (DirectionCosines, vtkMatrix4x4)
- [vtkGetObjectMacro](#) (Curve, vtkPolyData)
- [vtkGetVector3Macro](#) (ImagePositionPatient, double)
- [vtkGetVector6Macro](#) (ImageOrientationPatient, double)
- [vtkSetMacro](#) (LoadOverlays, int)
- [vtkSetMacro](#) (LoadIconImage, int)

- [vtkSetMacro](#) ([LossyFlag](#), int)
- [vtkSetMacro](#) ([ApplyLookupTable](#), int)
- [vtkTypeRevisionMacro](#) ([vtkGDCMImageReader2](#), [vtkMedicalImageReader2](#))

### Static Public Member Functions

- static [vtkGDCMImageReader2](#) \* [New](#) ()

### Protected Member Functions

- [vtkGDCMImageReader2](#) ()
- [~vtkGDCMImageReader2](#) ()
- void [FillMedicalImageInformation](#) (const [gdcmm::ImageReader](#) &reader)
- int [LoadSingleFile](#) (const char \*filename, char \*pointer, unsigned long &outlen)
- int [ProcessRequest](#) ([vtkInformation](#) \*request, [vtkInformationVector](#) \*\*inputVector, [vtkInformationVector](#) \*outputVector)
- int [RequestData](#) ([vtkInformation](#) \*request, [vtkInformationVector](#) \*\*inputVector, [vtkInformationVector](#) \*outputVector)
- int [RequestDataCompat](#) ()
- int [RequestInformation](#) ([vtkInformation](#) \*request, [vtkInformationVector](#) \*\*inputVector, [vtkInformationVector](#) \*outputVector)
- int [RequestInformationCompat](#) ()
- void [SetFilePattern](#) (const char \*)
- void [SetFilePrefix](#) (const char \*)
- [vtkGetStringMacro](#) ([FilePrefix](#))
- [vtkGetStringMacro](#) ([FilePattern](#))
- [vtkSetVector6Macro](#) ([ImageOrientationPatient](#), double)

### Protected Attributes

- int [ApplyInverseVideo](#)
- int [ApplyLookupTable](#)
- int [ApplyPlanarConfiguration](#)
- int [ApplyShiftScale](#)
- int [ApplyYBRToRGB](#)
- [vtkPolyData](#) \* [Curve](#)
- [vtkMatrix4x4](#) \* [DirectionCosines](#)
- int [ForceRescale](#)
- int [IconDataScalarType](#)
- int [IconImageDataExtent](#) [6]
- int [IconNumberOfScalarComponents](#)
- int [ImageFormat](#)
- double [ImageOrientationPatient](#) [6]
- double [ImagePositionPatient](#) [3]
- int [LoadIconImage](#)
- int [LoadOverlays](#)
- int [LossyFlag](#)
- int [NumberOfIconImages](#)
- int [NumberOfOverlays](#)
- int [PlanarConfiguration](#)
- double [Scale](#)
- double [Shift](#)

### 27.350.1 Detailed Description

Examples:

[Compute3DSpacing.cxx](#).

### 27.350.2 Constructor & Destructor Documentation

27.350.2.1 `vtkGDCMImageReader2::vtkGDCMImageReader2 ( )` [protected]

27.350.2.2 `vtkGDCMImageReader2::~~vtkGDCMImageReader2 ( )` [protected]

### 27.350.3 Member Function Documentation

27.350.3.1 `virtual int vtkGDCMImageReader2::CanReadFile ( const char * fname )` [virtual]

27.350.3.2 `void vtkGDCMImageReader2::FillMedicalImageInformation ( const gdcm::ImageReader & reader )`  
[protected]

27.350.3.3 `virtual const char* vtkGDCMImageReader2::GetDescriptiveName ( )` [inline],[virtual]

27.350.3.4 `virtual const char* vtkGDCMImageReader2::GetFileExtensions ( )` [inline],[virtual]

27.350.3.5 `vtkImageData* vtkGDCMImageReader2::GetIconImage ( )`

27.350.3.6 `vtkAlgorithmOutput* vtkGDCMImageReader2::GetIconImagePort ( )`

27.350.3.7 `vtkImageData* vtkGDCMImageReader2::GetOverlay ( int i )`

27.350.3.8 `vtkAlgorithmOutput* vtkGDCMImageReader2::GetOverlayPort ( int index )`

27.350.3.9 `int vtkGDCMImageReader2::LoadSingleFile ( const char * filename, char * pointer, unsigned long & outlen )`  
[protected]

27.350.3.10 `static vtkGDCMImageReader2* vtkGDCMImageReader2::New ( )` [static]

Examples:

[Compute3DSpacing.cxx](#).

27.350.3.11 `virtual void vtkGDCMImageReader2::PrintSelf ( ostream & os, vtkIndent indent )` [virtual]

27.350.3.12 `int vtkGDCMImageReader2::ProcessRequest ( vtkInformation * request, vtkInformationVector ** inputVector,  
vtkInformationVector * outputVector )` [protected]

27.350.3.13 `int vtkGDCMImageReader2::RequestData ( vtkInformation * request, vtkInformationVector ** inputVector,  
vtkInformationVector * outputVector )` [protected]

27.350.3.14 `int vtkGDCMImageReader2::RequestDataCompat ( )` [protected]

27.350.3.15 `int vtkGDCMImageReader2::RequestInformation ( vtkInformation * request, vtkInformationVector ** inputVector,  
vtkInformationVector * outputVector )` [protected]



- 27.350.3.16 int vtkGDCMImageReader2::RequestInformationCompat ( ) [protected]
- 27.350.3.17 virtual void vtkGDCMImageReader2::SetCurve ( vtkPolyData \* *pd* ) [virtual]
- 27.350.3.18 void vtkGDCMImageReader2::SetFilePattern ( const char \* ) [inline],[protected]
- 27.350.3.19 void vtkGDCMImageReader2::SetFilePrefix ( const char \* ) [inline],[protected]
- 27.350.3.20 virtual void vtkGDCMImageReader2::SetMedicalImageProperties ( vtkMedicalImageProperties \* *pd* ) [virtual]
- 27.350.3.21 vtkGDCMImageReader2::vtkBooleanMacro ( LoadOverlays , int )
- 27.350.3.22 vtkGDCMImageReader2::vtkBooleanMacro ( LoadIconImage , int )
- 27.350.3.23 vtkGDCMImageReader2::vtkBooleanMacro ( LossyFlag , int )
- 27.350.3.24 vtkGDCMImageReader2::vtkBooleanMacro ( ApplyLookupTable , int )
- 27.350.3.25 int vtkGDCMImageReader2::vtkBooleanMacro ( ApplyYBRToRGB , int )
- 27.350.3.26 vtkGDCMImageReader2::vtkGetMacro ( LoadOverlays , int )
- 27.350.3.27 vtkGDCMImageReader2::vtkGetMacro ( LoadIconImage , int )
- 27.350.3.28 vtkGDCMImageReader2::vtkGetMacro ( LossyFlag , int )
- 27.350.3.29 vtkGDCMImageReader2::vtkGetMacro ( NumberOfOverlays , int )
- 27.350.3.30 vtkGDCMImageReader2::vtkGetMacro ( NumberOfIconImages , int )
- 27.350.3.31 vtkGDCMImageReader2::vtkGetMacro ( ApplyLookupTable , int )
- 27.350.3.32 vtkGDCMImageReader2::vtkGetMacro ( ApplyYBRToRGB , int )
- 27.350.3.33 vtkGDCMImageReader2::vtkGetMacro ( ImageFormat , int )
- 27.350.3.34 vtkGDCMImageReader2::vtkGetMacro ( PlanarConfiguration , int )
- 27.350.3.35 vtkGDCMImageReader2::vtkGetMacro ( Shift , double )
- 27.350.3.36 vtkGDCMImageReader2::vtkGetMacro ( Scale , double )
- 27.350.3.37 vtkGDCMImageReader2::vtkGetObjectMacro ( DirectionCosines , vtkMatrix4x4 )
- 27.350.3.38 vtkGDCMImageReader2::vtkGetObjectMacro ( Curve , vtkPolyData )
- 27.350.3.39 vtkGDCMImageReader2::vtkGetStringMacro ( FilePrefix ) [protected]
- 27.350.3.40 vtkGDCMImageReader2::vtkGetStringMacro ( FilePattern ) [protected]
- 27.350.3.41 vtkGDCMImageReader2::vtkGetVector3Macro ( ImagePositionPatient , double )

- 27.350.3.42 `vtkGDCMImageReader2::vtkGetVector6Macro ( ImageOrientationPatient , double )`
- 27.350.3.43 `vtkGDCMImageReader2::vtkSetMacro ( LoadOverlays , int )`
- 27.350.3.44 `vtkGDCMImageReader2::vtkSetMacro ( LoadIconImage , int )`
- 27.350.3.45 `vtkGDCMImageReader2::vtkSetMacro ( LossyFlag , int )`
- 27.350.3.46 `vtkGDCMImageReader2::vtkSetMacro ( ApplyLookupTable , int )`
- 27.350.3.47 `vtkGDCMImageReader2::vtkSetVector6Macro ( ImageOrientationPatient , double )` [protected]
- 27.350.3.48 `vtkGDCMImageReader2::vtkTypeRevisionMacro ( vtkGDCMImageReader2 , vtkMedicalImageReader2 )`

#### 27.350.4 Member Data Documentation

- 27.350.4.1 `int vtkGDCMImageReader2::ApplyInverseVideo` [protected]
- 27.350.4.2 `int vtkGDCMImageReader2::ApplyLookupTable` [protected]
- 27.350.4.3 `int vtkGDCMImageReader2::ApplyPlanarConfiguration` [protected]
- 27.350.4.4 `int vtkGDCMImageReader2::ApplyShiftScale` [protected]
- 27.350.4.5 `int vtkGDCMImageReader2::ApplyYBRToRGB` [protected]
- 27.350.4.6 `vtkPolyData* vtkGDCMImageReader2::Curve` [protected]
- 27.350.4.7 `vtkMatrix4x4* vtkGDCMImageReader2::DirectionCosines` [protected]
- 27.350.4.8 `int vtkGDCMImageReader2::ForceRescale` [protected]
- 27.350.4.9 `int vtkGDCMImageReader2::IconDataScalarType` [protected]
- 27.350.4.10 `int vtkGDCMImageReader2::IconImageDataExtent[6]` [protected]
- 27.350.4.11 `int vtkGDCMImageReader2::IconNumberOfScalarComponents` [protected]
- 27.350.4.12 `int vtkGDCMImageReader2::ImageFormat` [protected]
- 27.350.4.13 `double vtkGDCMImageReader2::ImageOrientationPatient[6]` [protected]
- 27.350.4.14 `double vtkGDCMImageReader2::ImagePositionPatient[3]` [protected]
- 27.350.4.15 `int vtkGDCMImageReader2::LoadIconImage` [protected]
- 27.350.4.16 `int vtkGDCMImageReader2::LoadOverlays` [protected]
- 27.350.4.17 `int vtkGDCMImageReader2::LossyFlag` [protected]
- 27.350.4.18 `int vtkGDCMImageReader2::NumberOfIconImages` [protected]

27.350.4.19 int vtkGDCMImageReader2::NumberOfOverlays [protected]

27.350.4.20 int vtkGDCMImageReader2::PlanarConfiguration [protected]

27.350.4.21 double vtkGDCMImageReader2::Scale [protected]

27.350.4.22 double vtkGDCMImageReader2::Shift [protected]


The documentation for this class was generated from the following file:

- [vtkGDCMImageReader2.h](#)


## 27.351 vtkGDCMImageWriter Class Reference

```
#include <vtkGDCMImageWriter.h>
```

Inheritance diagram for vtkGDCMImageWriter:

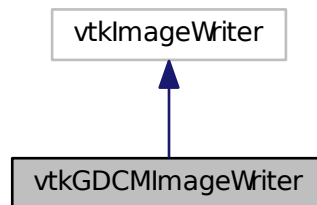

Collaboration diagram for vtkGDCMImageWriter:

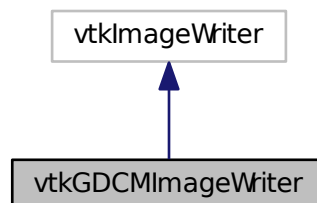

## Public Types

- enum [CompressionTypes](#) {  
[NO\\_COMPRESSION](#) = 0,  
[JPEG\\_COMPRESSION](#),  
[JPEG2000\\_COMPRESSION](#),  
[JPEGLS\\_COMPRESSION](#),  
[RLE\\_COMPRESSION](#) }

## Public Member Functions

- virtual const char \* [GetDescriptiveName](#) ()
- virtual const char \* [GetFileExtensions](#) ()
- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetDirectionCosines](#) (vtkMatrix4x4 \*matrix)
- virtual void [SetDirectionCosinesFromImageOrientationPatient](#) (const double dircos[6])
- virtual void [SetFileNames](#) (vtkStringArray \*)
- virtual void [SetMedicalImageProperties](#) (vtkMedicalImageProperties \*)
- [vtkBooleanMacro](#) (LossyFlag, int)
- [vtkBooleanMacro](#) (FileLowerLeft, int)
- [vtkGetMacro](#) (LossyFlag, int)
- [vtkGetMacro](#) (Shift, double)
- [vtkGetMacro](#) (Scale, double)
- [vtkGetMacro](#) (ImageFormat, int)
- [vtkGetMacro](#) (FileLowerLeft, int)
- [vtkGetMacro](#) (PlanarConfiguration, int)
- [vtkGetMacro](#) (CompressionType, int)
- [vtkGetObjectMacro](#) (MedicalImageProperties, vtkMedicalImageProperties)
- [vtkGetObjectMacro](#) (FileNames, vtkStringArray)
- [vtkGetObjectMacro](#) (DirectionCosines, vtkMatrix4x4)
- [vtkGetStringMacro](#) (StudyUID)
- [vtkGetStringMacro](#) (SeriesUID)
- [vtkSetMacro](#) (LossyFlag, int)
- [vtkSetMacro](#) (Shift, double)
- [vtkSetMacro](#) (Scale, double)
- [vtkSetMacro](#) (ImageFormat, int)
- [vtkSetMacro](#) (FileLowerLeft, int)
- [vtkSetMacro](#) (PlanarConfiguration, int)
- [vtkSetMacro](#) (CompressionType, int)
- [vtkSetStringMacro](#) (StudyUID)
- [vtkSetStringMacro](#) (SeriesUID)
- [vtkTypeRevisionMacro](#) (vtkGDCMImageWriter, vtkImageWriter)
- virtual void [Write](#) ()

## Static Public Member Functions

- static [vtkGDCMImageWriter](#) \* [New](#) ()

## Protected Member Functions

- [vtkGDCMImageWriter \(\)](#)
- [~vtkGDCMImageWriter \(\)](#)
- virtual char \* [GetFileName \(\)](#)
- int [WriteGDCMData](#) (vtkImageData \*data, int timeStep)
- void [WriteSlice](#) (vtkImageData \*data)

### 27.351.1 Detailed Description

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [CreateFakeRTDOSE.cxx](#), [gdcmorphoplanes.cxx](#), [HelloActiviz.cs](#), [HelloActiviz2.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld.java](#), [HelloVTKWorld2.cs](#), [MagnifyFile.cxx](#), and [RefCounting.cs](#).

### 27.351.2 Member Enumeration Documentation

#### 27.351.2.1 enum vtkGDCMImageWriter::CompressionTypes

Enumerator

***NO\_COMPRESSION***  
***JPEG\_COMPRESSION***  
***JPEG2000\_COMPRESSION***  
***JPEGLS\_COMPRESSION***  
***RLE\_COMPRESSION***

### 27.351.3 Constructor & Destructor Documentation

27.351.3.1 `vtkGDCMImageWriter::vtkGDCMImageWriter ( )` [protected]

27.351.3.2 `vtkGDCMImageWriter::~~vtkGDCMImageWriter ( )` [protected]

### 27.351.4 Member Function Documentation

27.351.4.1 `virtual const char* vtkGDCMImageWriter::GetDescriptiveName ( )` [inline],[virtual]

27.351.4.2 `virtual const char* vtkGDCMImageWriter::GetFileExtensions ( )` [inline],[virtual]

27.351.4.3 `virtual char* vtkGDCMImageWriter::GetFileName ( )` [protected],[virtual]

27.351.4.4 `static vtkGDCMImageWriter* vtkGDCMImageWriter::New ( )` [static]

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [CreateFakeRTDOSE.cxx](#), [gdcmorphoplanes.cxx](#), [HelloActiviz.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld2.cs](#), [MagnifyFile.cxx](#), and [RefCounting.cs](#).

27.351.4.5 virtual void vtkGDCMImageWriter::PrintSelf ( ostream & os, vtkIndent indent ) [virtual]

27.351.4.6 virtual void vtkGDCMImageWriter::SetDirectionCosines ( vtkMatrix4x4 \* matrix ) [virtual]

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [gdcmorphoplanes.cxx](#), [HelloActiviz2.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld.java](#), and [MagnifyFile.cxx](#).

27.351.4.7 virtual void vtkGDCMImageWriter::SetDirectionCosinesFromImageOrientationPatient ( const double dircos[6] ) [virtual]

27.351.4.8 virtual void vtkGDCMImageWriter::SetFileNames ( vtkStringArray \* ) [virtual]

Examples:

[ConvertMultiFrameToSingleFrame.cxx](#).

27.351.4.9 virtual void vtkGDCMImageWriter::SetMedicalImageProperties ( vtkMedicalImageProperties \* ) [virtual]

Examples:

[Convert16BitsTo8Bits.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [gdcmorphoplanes.cxx](#), [HelloActiviz.cs](#), [HelloActiviz2.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld.java](#), and [MagnifyFile.cxx](#).

27.351.4.10 vtkGDCMImageWriter::vtkBooleanMacro ( LossyFlag , int )

27.351.4.11 vtkGDCMImageWriter::vtkBooleanMacro ( FileLowerLeft , int )

27.351.4.12 vtkGDCMImageWriter::vtkGetMacro ( LossyFlag , int )

27.351.4.13 vtkGDCMImageWriter::vtkGetMacro ( Shift , double )

27.351.4.14 vtkGDCMImageWriter::vtkGetMacro ( Scale , double )

27.351.4.15 vtkGDCMImageWriter::vtkGetMacro ( ImageFormat , int )

27.351.4.16 vtkGDCMImageWriter::vtkGetMacro ( FileLowerLeft , int )

27.351.4.17 vtkGDCMImageWriter::vtkGetMacro ( PlanarConfiguration , int )

27.351.4.18 vtkGDCMImageWriter::vtkGetMacro ( CompressionType , int )

27.351.4.19 vtkGDCMImageWriter::vtkGetObjectMacro ( MedicalImageProperties , vtkMedicalImageProperties )

27.351.4.20 vtkGDCMImageWriter::vtkGetObjectMacro ( FileNames , vtkStringArray )

27.351.4.21 vtkGDCMImageWriter::vtkGetObjectMacro ( DirectionCosines , vtkMatrix4x4 )

27.351.4.22 vtkGDCMImageWriter::vtkGetStringMacro ( StudyUID )

- 27.351.4.23 `vtkGDCMImageWriter::vtkGetStringMacro ( SeriesUID )`
- 27.351.4.24 `vtkGDCMImageWriter::vtkSetMacro ( LossyFlag , int )`
- 27.351.4.25 `vtkGDCMImageWriter::vtkSetMacro ( Shift , double )`
- 27.351.4.26 `vtkGDCMImageWriter::vtkSetMacro ( Scale , double )`
- 27.351.4.27 `vtkGDCMImageWriter::vtkSetMacro ( ImageFormat , int )`
- 27.351.4.28 `vtkGDCMImageWriter::vtkSetMacro ( FileLowerLeft , int )`
- 27.351.4.29 `vtkGDCMImageWriter::vtkSetMacro ( PlanarConfiguration , int )`
- 27.351.4.30 `vtkGDCMImageWriter::vtkSetMacro ( CompressionType , int )`
- 27.351.4.31 `vtkGDCMImageWriter::vtkSetStringMacro ( StudyUID )`
- 27.351.4.32 `vtkGDCMImageWriter::vtkSetStringMacro ( SeriesUID )`
- 27.351.4.33 `vtkGDCMImageWriter::vtkTypeRevisionMacro ( vtkGDCMImageWriter , vtkImageWriter )`
- 27.351.4.34 `virtual void vtkGDCMImageWriter::Write ( ) [virtual]`

**Examples:**

[Convert16BitsTo8Bits.cxx](#), [ConvertMultiFrameToSingleFrame.cxx](#), [ConvertRGBToLuminance.cxx](#), [ConvertSingleBitTo8Bits.cxx](#), [CreateFakeRTDOSE.cxx](#), [gdcmorphoplanes.cxx](#), [HelloActiviz.cs](#), [HelloActiviz2.cs](#), [HelloVTKWorld.cs](#), [HelloVTKWorld.java](#), [HelloVTKWorld2.cs](#), and [MagnifyFile.cxx](#).

- 27.351.4.35 `int vtkGDCMImageWriter::WriteGDCMData ( vtkImageData * data, int timeStep ) [protected]`
- 27.351.4.36 `void vtkGDCMImageWriter::WriteSlice ( vtkImageData * data ) [protected]`


The documentation for this class was generated from the following file:

- [vtkGDCMImageWriter.h](#)


## 27.352 vtkGDCMMedicalImageProperties Class Reference

```
#include <vtkGDCMMedicalImageProperties.h>
```

Inheritance diagram for `vtkGDCMMedicalImageProperties`:

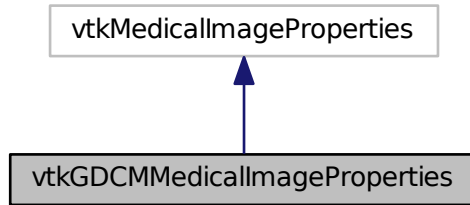

Collaboration diagram for `vtkGDCMMedicalImageProperties`:

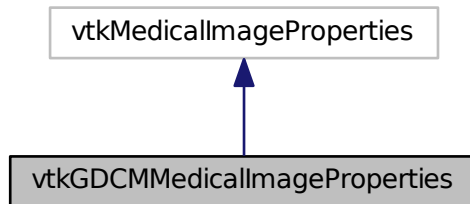

### Public Member Functions

- virtual void [Clear](#) ()
- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkTypeRevisionMacro](#) ([vtkGDCMMedicalImageProperties](#), vtkMedicalImageProperties)

### Static Public Member Functions

- static [vtkGDCMMedicalImageProperties](#) \* [New](#) ()

### Protected Member Functions

- [vtkGDCMMedicalImageProperties](#) ()
- [~vtkGDCMMedicalImageProperties](#) ()
- [gdcmm::File](#) const & [GetFile](#) (unsigned int t)
- void [PushBackFile](#) ([gdcmm::File](#) const &f)



## Friends

- class [vtkGDCMImageReader](#)
- class [vtkGDCMImageReader2](#)
- class [vtkGDCMImageWriter](#)

## 27.352.1 Constructor & Destructor Documentation

27.352.1.1 `vtkGDCMMedicalImageProperties::vtkGDCMMedicalImageProperties ( )` [protected]

27.352.1.2 `vtkGDCMMedicalImageProperties::~~vtkGDCMMedicalImageProperties ( )` [protected]

## 27.352.2 Member Function Documentation

27.352.2.1 `virtual void vtkGDCMMedicalImageProperties::Clear ( )` [virtual]

27.352.2.2 `gdcmm::File const& vtkGDCMMedicalImageProperties::GetFile ( unsigned int t )` [protected]

27.352.2.3 `static vtkGDCMMedicalImageProperties* vtkGDCMMedicalImageProperties::New ( )` [static]

27.352.2.4 `void vtkGDCMMedicalImageProperties::PrintSelf ( ostream & os, vtkIndent indent )`

27.352.2.5 `void vtkGDCMMedicalImageProperties::PushBackFile ( gdcmm::File const & f )` [protected]

27.352.2.6 `vtkGDCMMedicalImageProperties::vtkTypeRevisionMacro ( vtkGDCMMedicalImageProperties ,  
vtkMedicalImageProperties )`

## 27.352.3 Friends And Related Function Documentation

27.352.3.1 `friend class vtkGDCMImageReader` [friend]

27.352.3.2 `friend class vtkGDCMImageReader2` [friend]

27.352.3.3 `friend class vtkGDCMImageWriter` [friend]


The documentation for this class was generated from the following file:

- [vtkGDCMMedicalImageProperties.h](#)


## 27.353 vtkGDCMPolyDataReader Class Reference

```
#include <vtkGDCMPolyDataReader.h>
```

Inheritance diagram for vtkGDCMPolyDataReader:

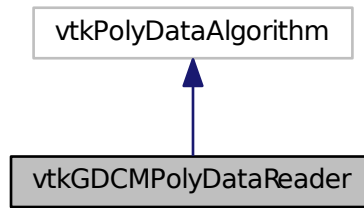

Collaboration diagram for vtkGDCMPolyDataReader:

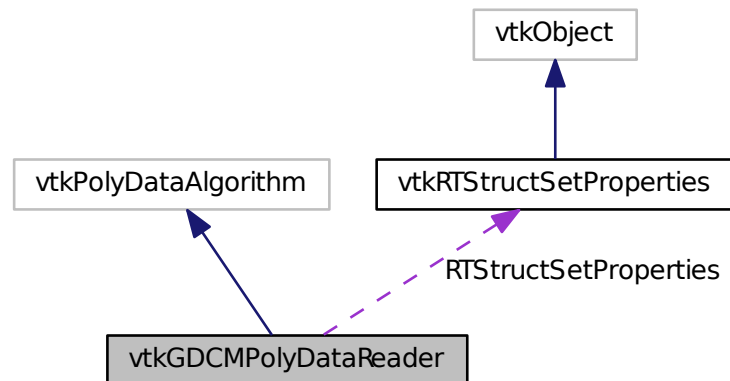

## Public Member Functions

- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkGetObjectMacro](#) ([MedicalImageProperties](#), vtkMedicalImageProperties)
- [vtkGetObjectMacro](#) ([RTStructSetProperties](#), [vtkRTStructSetProperties](#))
- [vtkGetStringMacro](#) ([FileName](#))
- [vtkSetStringMacro](#) ([FileName](#))
- [vtkTypeRevisionMacro](#) ([vtkGDCMPolyDataReader](#), vtkPolyDataAlgorithm)

## Static Public Member Functions

- static [vtkGDCMPolyDataReader](#) \* [New](#) ()

## Protected Member Functions

- [vtkGDCMPolyDataReader \(\)](#)
- [~vtkGDCMPolyDataReader \(\)](#)
- void [FillMedicalImageInformation](#) (const [gdcmm::Reader](#) &reader)
- int [RequestData](#) (vtkInformation \*, vtkInformationVector \*\*, vtkInformationVector \*)
- int [RequestData\\_HemodynamicWaveformStorage](#) ([gdcmm::Reader](#) const &reader, vtkInformationVector \*outputVector)
- int [RequestData\\_RTStructureSetStorage](#) ([gdcmm::Reader](#) const &reader, vtkInformationVector \*outputVector)
- int [RequestInformation](#) (vtkInformation \*vtkNotUsed(request), vtkInformationVector \*\*vtkNotUsed(inputVector), vtkInformationVector \*outputVector)
- int [RequestInformation\\_HemodynamicWaveformStorage](#) ([gdcmm::Reader](#) const &reader)
- int [RequestInformation\\_RTStructureSetStorage](#) ([gdcmm::Reader](#) const &reader)

## Protected Attributes

- char \* [FileName](#)
- vtkMedicalImageProperties \* [MedicalImageProperties](#)
- [vtkRTStructSetProperties](#) \* [RTStructSetProperties](#)

### 27.353.1 Detailed Description

Examples:

[gdcmscene.cxx](#), [GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

### 27.353.2 Constructor & Destructor Documentation

27.353.2.1 [vtkGDCMPolyDataReader::vtkGDCMPolyDataReader \( \)](#) [protected]

27.353.2.2 [vtkGDCMPolyDataReader::~~vtkGDCMPolyDataReader \( \)](#) [protected]

### 27.353.3 Member Function Documentation

27.353.3.1 void [vtkGDCMPolyDataReader::FillMedicalImageInformation \( const \[gdcmm::Reader\]\(#\) & reader \)](#) [protected]

27.353.3.2 static [vtkGDCMPolyDataReader\\*](#) [vtkGDCMPolyDataReader::New \( \)](#) [static]

Examples:

[gdcmscene.cxx](#), [GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

27.353.3.3 virtual void [vtkGDCMPolyDataReader::PrintSelf \( ostream & os, vtkIndent indent \)](#) [virtual]

27.353.3.4 int [vtkGDCMPolyDataReader::RequestData \( vtkInformation \\*, vtkInformationVector \\*\\*, vtkInformationVector \\* \)](#) [protected]

27.353.3.5 int [vtkGDCMPolyDataReader::RequestData\\_HemodynamicWaveformStorage \( \[gdcmm::Reader\]\(#\) const & reader, vtkInformationVector \\* outputVector \)](#) [protected]

- 27.353.3.6 `int vtkGDCMPolyDataReader::RequestData_RTStructureSetStorage ( gdcm::Reader const & reader, vtkInformationVector * outputVector )` [protected]
- 27.353.3.7 `int vtkGDCMPolyDataReader::RequestInformation ( vtkInformation * vtkNotUsedrequest, vtkInformationVector ** vtkNotUsedinputVector, vtkInformationVector * outputVector )` [protected]
- 27.353.3.8 `int vtkGDCMPolyDataReader::RequestInformation_HemodynamicWaveformStorage ( gdcm::Reader const & reader )` [protected]
- 27.353.3.9 `int vtkGDCMPolyDataReader::RequestInformation_RTStructureSetStorage ( gdcm::Reader const & reader )` [protected]
- 27.353.3.10 `vtkGDCMPolyDataReader::vtkGetObjectMacro ( MedicalImageProperties , vtkMedicalImageProperties )`
- 27.353.3.11 `vtkGDCMPolyDataReader::vtkGetObjectMacro ( RTStructSetProperties , vtkRTStructSetProperties )`
- 27.353.3.12 `vtkGDCMPolyDataReader::vtkGetStringMacro ( FileName )`
- 27.353.3.13 `vtkGDCMPolyDataReader::vtkSetStringMacro ( FileName )`
- 27.353.3.14 `vtkGDCMPolyDataReader::vtkTypeRevisionMacro ( vtkGDCMPolyDataReader , vtkPolyDataAlgorithm )`

## 27.353.4 Member Data Documentation

- 27.353.4.1 `char* vtkGDCMPolyDataReader::FileName` [protected]
- 27.353.4.2 `vtkMedicalImageProperties* vtkGDCMPolyDataReader::MedicalImageProperties` [protected]
- 27.353.4.3 `vtkRTStructSetProperties* vtkGDCMPolyDataReader::RTStructSetProperties` [protected]


The documentation for this class was generated from the following file:

- [vtkGDCMPolyDataReader.h](#)


## 27.354 vtkGDCMPolyDataWriter Class Reference

```
#include <vtkGDCMPolyDataWriter.h>
```

Inheritance diagram for vtkGDCMPolyDataWriter:

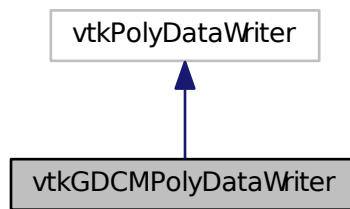

Collaboration diagram for vtkGDCMPolyDataWriter:

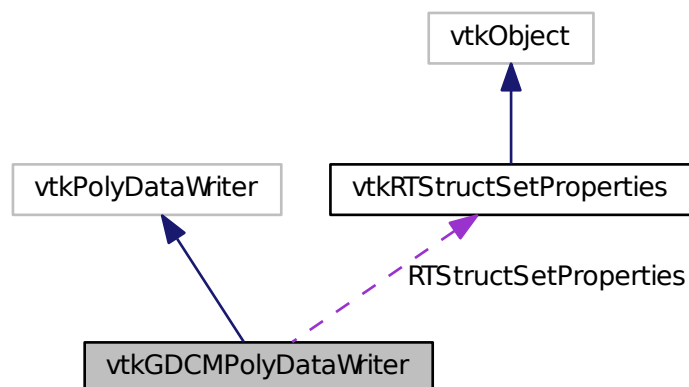

## Public Member Functions

- void [InitializeRTStructSet](#) (vtkStdString inDirectory, vtkStdString inStructLabel, vtkStdString inStructName, vtkStringArray \*inROINames, vtkStringArray \*inROIAlgorithmName, vtkStringArray \*inROIType)
- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetMedicalImageProperties](#) (vtkMedicalImageProperties \*pd)
- void [SetNumberOfInputPorts](#) (int n)
- virtual void [SetRTStructSetProperties](#) (vtkRTStructSetProperties \*pd)
- [vtkTypeRevisionMacro](#) (vtkGDCMPolyDataWriter, vtkPolyDataWriter)

## Static Public Member Functions

- static [vtkGDCMPolyDataWriter \\* New](#) ()

## Protected Member Functions

- [vtkGDCMPolyDataWriter](#) ()
- [~vtkGDCMPolyDataWriter](#) ()
- void [WriteData](#) ()
- void [WriteRTSTRUCTData](#) ([gdcmm::File](#) &file, int num)
- void [WriteRTSTRUCTInfo](#) ([gdcmm::File](#) &file)

## Protected Attributes

- [vtkMedicalImageProperties](#) \* [MedicalImageProperties](#)
- [vtkRTStructSetProperties](#) \* [RTStructSetProperties](#)

### 27.354.1 Detailed Description

Examples:

[GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

### 27.354.2 Constructor & Destructor Documentation

27.354.2.1 `vtkGDCMPolyDataWriter::vtkGDCMPolyDataWriter ( )` [protected]

27.354.2.2 `vtkGDCMPolyDataWriter::~~vtkGDCMPolyDataWriter ( )` [protected]

### 27.354.3 Member Function Documentation

27.354.3.1 `void vtkGDCMPolyDataWriter::InitializeRTStructSet ( vtkStdString inDirectory, vtkStdString inStructLabel, vtkStdString inStructName, vtkStringArray * inROINames, vtkStringArray * inROIAlgorithmName, vtkStringArray * inROIType )`

Examples:

[GenerateRTSTRUCT.cxx](#).

27.354.3.2 `static vtkGDCMPolyDataWriter* vtkGDCMPolyDataWriter::New ( )` [static]

Examples:

[GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

27.354.3.3 `virtual void vtkGDCMPolyDataWriter::PrintSelf ( ostream & os, vtkIndent indent )` [virtual]

27.354.3.4 `virtual void vtkGDCMPolyDataWriter::SetMedicalImageProperties ( vtkMedicalImageProperties * pd )` [virtual]

Examples:

[GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

27.354.3.5 void vtkGDCMPolyDataWriter::SetNumberOfInputPorts ( int *n* )

Examples:

[GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

27.354.3.6 virtual void vtkGDCMPolyDataWriter::SetRTStructSetProperties ( vtkRTStructSetProperties \* *pd* ) [virtual]

Examples:

[GenerateRTSTRUCT.cxx](#), and [rtstructapp.cxx](#).

27.354.3.7 vtkGDCMPolyDataWriter::vtkTypeRevisionMacro ( vtkGDCMPolyDataWriter , vtkPolyDataWriter )

27.354.3.8 void vtkGDCMPolyDataWriter::WriteData ( ) [protected]

27.354.3.9 void vtkGDCMPolyDataWriter::WriteRTSTRUCTData ( gdcm::File & *file*, int *num* ) [protected]

27.354.3.10 void vtkGDCMPolyDataWriter::WriteRTSTRUCTInfo ( gdcm::File & *file* ) [protected]

## 27.354.4 Member Data Documentation

27.354.4.1 vtkMedicalImageProperties\* vtkGDCMPolyDataWriter::MedicalImageProperties [protected]

27.354.4.2 vtkRTStructSetProperties\* vtkGDCMPolyDataWriter::RTStructSetProperties [protected]


The documentation for this class was generated from the following file:

- [vtkGDCMPolyDataWriter.h](#)


## 27.355 vtkGDCMTesting Class Reference

```
#include <vtkGDCMTesting.h>
```

Inheritance diagram for vtkGDCMTesting:

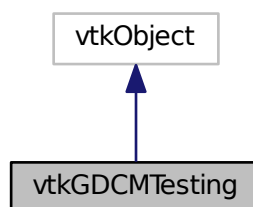

Collaboration diagram for vtkGDCMTesting:

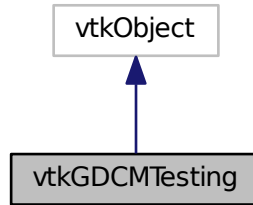

## Public Types

- typedef const char \*const (\* [MD5MetalImagesType](#))[3]

## Public Member Functions

- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkTypeRevisionMacro](#) (vtkGDCMTesting, vtkObject)

## Static Public Member Functions

- static const char \* [GetGDCMDataRoot](#) ()
- static const char \*const \* [GetMD5MetalImage](#) (unsigned int file)
- static const char \* [GetMHDMD5FromFile](#) (const char \*filepath)
- static unsigned int [GetNumberOfMD5MetalImages](#) ()
- static const char \* [GetRAWMD5FromFile](#) (const char \*filepath)
- static const char \* [GetVTKDataRoot](#) ()
- static [vtkGDCMTesting](#) \* [New](#) ()

## Protected Member Functions

- [vtkGDCMTesting](#) ()
- [~vtkGDCMTesting](#) ()

## 27.355.1 Detailed Description

Examples:

[HelloActiviz5.cs](#), [HelloVTKWorld2.cs](#), [MetalImageMD5Activiz.cs](#), [ReadSeriesIntoVTK.java](#), and [RefCounting.cs](#).



## 27.355.2 Member Typedef Documentation

27.355.2.1 `typedef const char* const(* vtkGDCMTesting::MD5MetalmagesType)[3]`

## 27.355.3 Constructor & Destructor Documentation

27.355.3.1 `vtkGDCMTesting::vtkGDCMTesting ( )` [protected]

27.355.3.2 `vtkGDCMTesting::~~vtkGDCMTesting ( )` [protected]

## 27.355.4 Member Function Documentation

27.355.4.1 `static const char* vtkGDCMTesting::GetGDCMDataRoot ( )` [static]

Examples:

[HelloActiviz5.cs](#), and [ReadSeriesIntoVTK.java](#).

27.355.4.2 `static const char* const* vtkGDCMTesting::GetMD5Metalmage ( unsigned int file )` [static]

27.355.4.3 `static const char* vtkGDCMTesting::GetMHDM5FromFile ( const char * filepath )` [static]

Examples:

[MetalmageMD5Activiz.cs](#).

27.355.4.4 `static unsigned int vtkGDCMTesting::GetNumberOfMD5Metalmages ( )` [static]

27.355.4.5 `static const char* vtkGDCMTesting::GetRAWMD5FromFile ( const char * filepath )` [static]

Examples:

[MetalmageMD5Activiz.cs](#).

27.355.4.6 `static const char* vtkGDCMTesting::GetVTKDataRoot ( )` [static]

Examples:

[HelloActiviz5.cs](#), and [HelloVTKWorld2.cs](#).

27.355.4.7 `static vtkGDCMTesting* vtkGDCMTesting::New ( )` [static]

Examples:

[RefCounting.cs](#).

27.355.4.8 void vtkGDCMTesting::PrintSelf ( ostream & *os*, vtkIndent *indent* )

27.355.4.9 vtkGDCMTesting::vtkTypeRevisionMacro ( vtkGDCMTesting , vtkObject )


The documentation for this class was generated from the following file:

- [vtkGDCMTesting.h](#)


## 27.356 vtkGDCMThreadedImageReader Class Reference

```
#include <vtkGDCMThreadedImageReader.h>
```

Inheritance diagram for vtkGDCMThreadedImageReader:

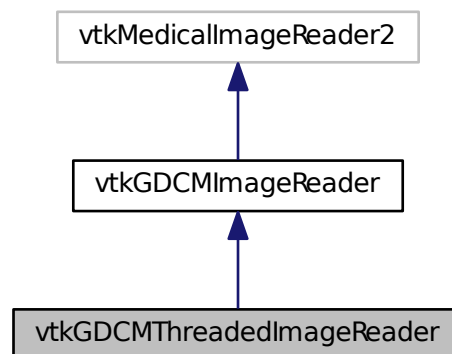

Collaboration diagram for vtkGDCMThreadedImageReader:

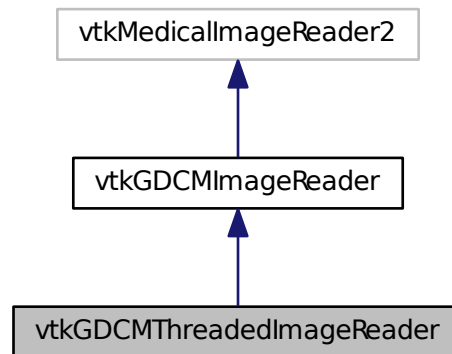

### Public Member Functions

- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkBooleanMacro](#) (UseShiftScale, int)
- [vtkGetMacro](#) (UseShiftScale, int)
- [vtkSetMacro](#) (Shift, double)
- [vtkSetMacro](#) (Scale, double)
- [vtkSetMacro](#) (UseShiftScale, int)
- [vtkTypeRevisionMacro](#) (vtkGDCMThreadedImageReader, vtkGDCMImageReader)

### Static Public Member Functions

- static [vtkGDCMThreadedImageReader \\* New](#) ()

### Protected Member Functions

- [vtkGDCMThreadedImageReader](#) ()
- [~vtkGDCMThreadedImageReader](#) ()
- void [ExecuteData](#) (vtkDataObject \*out)
- void [ExecuteInformation](#) ()
- void [ReadFiles](#) (unsigned int nfiles, const char \*filenames[])
- void [RequestDataCompat](#) ()

### Additional Inherited Members

#### 27.356.1 Constructor & Destructor Documentation

27.356.1.1 [vtkGDCMThreadedImageReader::vtkGDCMThreadedImageReader](#) ( ) `[protected]`

27.356.1.2 `vtkGDCMThreadedImageReader::~~vtkGDCMThreadedImageReader ( )` [protected]

## 27.356.2 Member Function Documentation

27.356.2.1 `void vtkGDCMThreadedImageReader::ExecuteData ( vtkDataObject * out )` [protected]

27.356.2.2 `void vtkGDCMThreadedImageReader::ExecuteInformation ( )` [protected]

27.356.2.3 `static vtkGDCMThreadedImageReader* vtkGDCMThreadedImageReader::New ( )` [static]

27.356.2.4 `virtual void vtkGDCMThreadedImageReader::PrintSelf ( ostream & os, vtkIndent indent )` [virtual]

Reimplemented from [vtkGDCMImageReader](#).

27.356.2.5 `void vtkGDCMThreadedImageReader::ReadFiles ( unsigned int nfiles, const char * filenames[ ] )` [protected]

27.356.2.6 `void vtkGDCMThreadedImageReader::RequestDataCompat ( )` [protected]

27.356.2.7 `vtkGDCMThreadedImageReader::vtkBooleanMacro ( UseShiftScale , int )`

27.356.2.8 `vtkGDCMThreadedImageReader::vtkGetMacro ( UseShiftScale , int )`

27.356.2.9 `vtkGDCMThreadedImageReader::vtkSetMacro ( Shift , double )`

27.356.2.10 `vtkGDCMThreadedImageReader::vtkSetMacro ( Scale , double )`

27.356.2.11 `vtkGDCMThreadedImageReader::vtkSetMacro ( UseShiftScale , int )`

27.356.2.12 `vtkGDCMThreadedImageReader::vtkTypeRevisionMacro ( vtkGDCMThreadedImageReader ,  
vtkGDCMImageReader )`


The documentation for this class was generated from the following file:

- [vtkGDCMThreadedImageReader.h](#)


## 27.357 vtkGDCMThreadedImageReader2 Class Reference

```
#include <vtkGDCMThreadedImageReader2.h>
```

Inheritance diagram for vtkGDCMThreadedImageReader2:

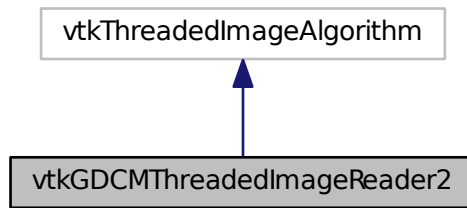

Collaboration diagram for vtkGDCMThreadedImageReader2:

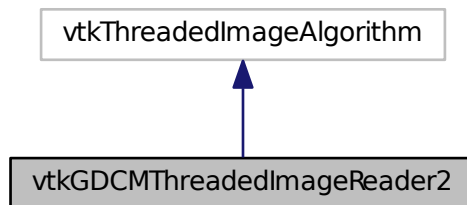

## Public Member Functions

- virtual const char \* [GetFileName](#) (int i=0)
- virtual void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetFileName](#) (const char \*filename)
- virtual void [SetFileNames](#) (vtkStringArray \*)
- int [SplitExtent](#) (int splitExt[6], int startExt[6], int num, int total)
- [vtkBooleanMacro](#) (FileLowerLeft, int)
- [vtkBooleanMacro](#) (LoadOverlays, int)
- [vtkBooleanMacro](#) (UseShiftScale, int)
- [vtkGetMacro](#) (FileLowerLeft, int)
- [vtkGetMacro](#) (NumberOfOverlays, int)
- [vtkGetMacro](#) (DataScalarType, int)
- [vtkGetMacro](#) (NumberOfScalarComponents, int)
- [vtkGetMacro](#) (LoadOverlays, int)
- [vtkGetMacro](#) (Shift, double)
- [vtkGetMacro](#) (Scale, double)
- [vtkGetMacro](#) (UseShiftScale, int)

- [vtkGetObjectMacro](#) (FileNames, vtkStringArray)
- [vtkGetVector3Macro](#) (DataOrigin, double)
- [vtkGetVector3Macro](#) (DataSpacing, double)
- [vtkGetVector6Macro](#) (DataExtent, int)
- [vtkSetMacro](#) (FileLowerLeft, int)
- [vtkSetMacro](#) (DataScalarType, int)
- [vtkSetMacro](#) (NumberOfScalarComponents, int)
- [vtkSetMacro](#) (LoadOverlays, int)
- [vtkSetMacro](#) (Shift, double)
- [vtkSetMacro](#) (Scale, double)
- [vtkSetMacro](#) (UseShiftScale, int)
- [vtkSetVector3Macro](#) (DataOrigin, double)
- [vtkSetVector3Macro](#) (DataSpacing, double)
- [vtkSetVector6Macro](#) (DataExtent, int)
- [vtkTypeRevisionMacro](#) ([vtkGDCMThreadedImageReader2](#), vtkThreadedImageAlgorithm)

### Static Public Member Functions

- static [vtkGDCMThreadedImageReader2](#) \* [New](#) ()

### Protected Member Functions

- [vtkGDCMThreadedImageReader2](#) ()
- [~vtkGDCMThreadedImageReader2](#) ()
- int [RequestInformation](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector)
- void [ThreadedRequestData](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector, vtkImageData \*\*\*inData, vtkImageData \*\*outData, int outExt[6], int id)

## 27.357.1 Constructor & Destructor Documentation

27.357.1.1 [vtkGDCMThreadedImageReader2::vtkGDCMThreadedImageReader2](#) ( ) [protected]

27.357.1.2 [vtkGDCMThreadedImageReader2::~~vtkGDCMThreadedImageReader2](#) ( ) [protected]

## 27.357.2 Member Function Documentation

27.357.2.1 virtual const char\* [vtkGDCMThreadedImageReader2::GetFileName](#) ( int *i* = 0 ) [virtual]

27.357.2.2 static [vtkGDCMThreadedImageReader2](#)\* [vtkGDCMThreadedImageReader2::New](#) ( ) [static]

27.357.2.3 virtual void [vtkGDCMThreadedImageReader2::PrintSelf](#) ( ostream & *os*, vtkIndent *indent* ) [virtual]

27.357.2.4 int [vtkGDCMThreadedImageReader2::RequestInformation](#) ( vtkInformation \* *request*, vtkInformationVector \*\* *inputVector*, vtkInformationVector \* *outputVector* ) [protected]

27.357.2.5 virtual void [vtkGDCMThreadedImageReader2::SetFileName](#) ( const char \* *filename* ) [virtual]

27.357.2.6 virtual void [vtkGDCMThreadedImageReader2::SetFileNames](#) ( vtkStringArray \* ) [virtual]

- 27.357.2.7 `int vtkGDCMThreadedImageReader2::SplitExtent ( int splitExt[6], int startExt[6], int num, int total )`
- 27.357.2.8 `void vtkGDCMThreadedImageReader2::ThreadedRequestData ( vtkInformation * request, vtkInformationVector ** inputVector, vtkInformationVector * outputVector, vtkImageData *** inData, vtkImageData ** outData, int outExt[6], int id ) [protected]`
- 27.357.2.9 `vtkGDCMThreadedImageReader2::vtkBooleanMacro ( FileLowerLeft , int )`
- 27.357.2.10 `vtkGDCMThreadedImageReader2::vtkBooleanMacro ( LoadOverlays , int )`
- 27.357.2.11 `vtkGDCMThreadedImageReader2::vtkBooleanMacro ( UseShiftScale , int )`
- 27.357.2.12 `vtkGDCMThreadedImageReader2::vtkGetMacro ( FileLowerLeft , int )`
- 27.357.2.13 `vtkGDCMThreadedImageReader2::vtkGetMacro ( NumberOfOverlays , int )`
- 27.357.2.14 `vtkGDCMThreadedImageReader2::vtkGetMacro ( DataScalarType , int )`
- 27.357.2.15 `vtkGDCMThreadedImageReader2::vtkGetMacro ( NumberOfScalarComponents , int )`
- 27.357.2.16 `vtkGDCMThreadedImageReader2::vtkGetMacro ( LoadOverlays , int )`
- 27.357.2.17 `vtkGDCMThreadedImageReader2::vtkGetMacro ( Shift , double )`
- 27.357.2.18 `vtkGDCMThreadedImageReader2::vtkGetMacro ( Scale , double )`
- 27.357.2.19 `vtkGDCMThreadedImageReader2::vtkGetMacro ( UseShiftScale , int )`
- 27.357.2.20 `vtkGDCMThreadedImageReader2::vtkGetObjectMacro ( FileNames , vtkStringArray )`
- 27.357.2.21 `vtkGDCMThreadedImageReader2::vtkGetVector3Macro ( DataOrigin , double )`
- 27.357.2.22 `vtkGDCMThreadedImageReader2::vtkGetVector3Macro ( DataSpacing , double )`
- 27.357.2.23 `vtkGDCMThreadedImageReader2::vtkGetVector6Macro ( DataExtent , int )`
- 27.357.2.24 `vtkGDCMThreadedImageReader2::vtkSetMacro ( FileLowerLeft , int )`
- 27.357.2.25 `vtkGDCMThreadedImageReader2::vtkSetMacro ( DataScalarType , int )`
- 27.357.2.26 `vtkGDCMThreadedImageReader2::vtkSetMacro ( NumberOfScalarComponents , int )`
- 27.357.2.27 `vtkGDCMThreadedImageReader2::vtkSetMacro ( LoadOverlays , int )`
- 27.357.2.28 `vtkGDCMThreadedImageReader2::vtkSetMacro ( Shift , double )`
- 27.357.2.29 `vtkGDCMThreadedImageReader2::vtkSetMacro ( Scale , double )`
- 27.357.2.30 `vtkGDCMThreadedImageReader2::vtkSetMacro ( UseShiftScale , int )`
- 27.357.2.31 `vtkGDCMThreadedImageReader2::vtkSetVector3Macro ( DataOrigin , double )`

27.357.2.32 `vtkGDCMThreadedImageReader2::vtkSetVector3Macro ( DataSpacing , double )`

27.357.2.33 `vtkGDCMThreadedImageReader2::vtkSetVector6Macro ( DataExtent , int )`

27.357.2.34 `vtkGDCMThreadedImageReader2::vtkTypeRevisionMacro ( vtkGDCMThreadedImageReader2 ,  
vtkThreadedImageAlgorithm )`


The documentation for this class was generated from the following file:

- [vtkGDCMThreadedImageReader2.h](#)

## 27.358 vtkImageColorViewer Class Reference


```
#include <vtkImageColorViewer.h>
```

Inheritance diagram for `vtkImageColorViewer`:

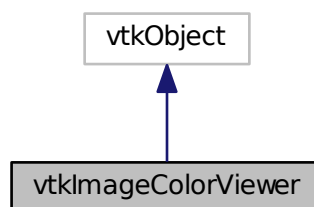



Collaboration diagram for vtkImageColorViewer:

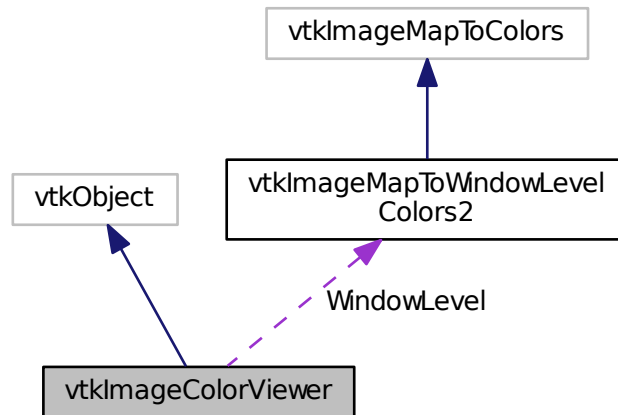

## Public Types

- enum {  
[SLICE\\_ORIENTATION\\_YZ](#) = 0,  
[SLICE\\_ORIENTATION\\_XZ](#) = 1,  
[SLICE\\_ORIENTATION\\_XY](#) = 2 }

## Public Member Functions

- virtual void [AddInput](#) (vtkImageData \*input)
- virtual void [AddInputConnection](#) (vtkAlgorithmOutput \*input)
- virtual double [GetColorLevel](#) ()
- virtual double [GetColorWindow](#) ()
- virtual vtkImageData \* [GetInput](#) ()
- virtual int [GetOffScreenRendering](#) ()
- double [GetOverlayVisibility](#) ()
- virtual int \* [GetPosition](#) ()
- virtual int \* [GetSize](#) ()
- virtual int [GetSliceMax](#) ()
- virtual int [GetSliceMin](#) ()
- virtual void [GetSliceRange](#) (int range[2])
- virtual void [GetSliceRange](#) (int &min, int &max)
- virtual int \* [GetSliceRange](#) ()
- virtual const char \* [GetWindowName](#) ()
- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [Render](#) (void)
- virtual void [SetColorLevel](#) (double s)
- virtual void [SetColorWindow](#) (double s)

- virtual void [SetDisplayId](#) (void \*a)
- virtual void [SetInput](#) (vtkImageData \*in)
- virtual void [SetInputConnection](#) (vtkAlgorithmOutput \*input)
- virtual void [SetOffScreenRendering](#) (int)
- void [SetOverlayVisibility](#) (double vis)
- virtual void [SetParentId](#) (void \*a)
- virtual void [SetPosition](#) (int a, int b)
- virtual void [SetPosition](#) (int a[2])
- virtual void [SetRenderer](#) (vtkRenderer \*arg)
- virtual void [SetRenderWindow](#) (vtkRenderWindow \*arg)
- virtual void [SetSize](#) (int a, int b)
- virtual void [SetSize](#) (int a[2])
- virtual void [SetSlice](#) (int s)
- virtual void [SetSliceOrientation](#) (int orientation)
- virtual void [SetSliceOrientationToXY](#) ()
- virtual void [SetSliceOrientationToXZ](#) ()
- virtual void [SetSliceOrientationToYZ](#) ()
- virtual void [SetupInteractor](#) (vtkRenderWindowInteractor \*)
- virtual void [SetWindowId](#) (void \*a)
- virtual void [UpdateDisplayExtent](#) ()
- [VTK\\_LEGACY](#) (int GetWholeZMin())
- [VTK\\_LEGACY](#) (int GetWholeZMax())
- [VTK\\_LEGACY](#) (int GetZSlice())
- [VTK\\_LEGACY](#) (void SetZSlice(int))
- [vtkBooleanMacro](#) (OffScreenRendering, int)
- [vtkGetMacro](#) (SliceOrientation, int)
- [vtkGetMacro](#) (Slice, int)
- [vtkGetObjectMacro](#) (RenderWindow, vtkRenderWindow)
- [vtkGetObjectMacro](#) (Renderer, vtkRenderer)
- [vtkGetObjectMacro](#) (ImageActor, vtkImageActor)
- [vtkGetObjectMacro](#) (WindowLevel, vtkImageMapToWindowLevelColors2)
- [vtkGetObjectMacro](#) (InteractorStyle, vtkInteractorStyleImage)
- [vtkTypeRevisionMacro](#) (vtkImageColorViewer, vtkObject)

## Static Public Member Functions

- static [vtkImageColorViewer](#) \* [New](#) ()

## Protected Member Functions

- [vtkImageColorViewer](#) ()
- [~vtkImageColorViewer](#) ()
- virtual void [InstallPipeline](#) ()
- virtual void [UnInstallPipeline](#) ()
- virtual void [UpdateOrientation](#) ()

## Protected Attributes

- int [FirstRender](#)
- vtkImageActor \* [ImageActor](#)
- vtkRenderWindowInteractor \* [Interactor](#)
- vtkInteractorStyleImage \* [InteractorStyle](#)
- vtkImageActor \* [OverlayImageActor](#)
- vtkRenderer \* [Renderer](#)
- vtkRenderWindow \* [RenderWindow](#)
- int [Slice](#)
- int [SliceOrientation](#)
- [vtkImageMapToWindowLevelColors2](#) \* [WindowLevel](#)

## Friends

- class [vtkImageColorViewerCallback](#)

### 27.358.1 Detailed Description

Examples:

[gdcmrtionplan.cxx](#), and [gdcmrtplan.cxx](#).

### 27.358.2 Member Enumeration Documentation

27.358.2.1 anonymous enum

Enumerator

***SLICE\_ORIENTATION\_YZ***  
***SLICE\_ORIENTATION\_XZ***  
***SLICE\_ORIENTATION\_XY***

### 27.358.3 Constructor & Destructor Documentation

27.358.3.1 `vtkImageColorViewer::vtkImageColorViewer ( )` [protected]

27.358.3.2 `vtkImageColorViewer::~~vtkImageColorViewer ( )` [protected]

### 27.358.4 Member Function Documentation

27.358.4.1 `virtual void vtkImageColorViewer::AddInput ( vtkImageData * input )` [virtual]

27.358.4.2 `virtual void vtkImageColorViewer::AddInputConnection ( vtkAlgorithmOutput * input )` [virtual]

27.358.4.3 `virtual double vtkImageColorViewer::GetColorLevel ( )` [virtual]

27.358.4.4 `virtual double vtkImageColorViewer::GetColorWindow ( )` [virtual]

27.358.4.5 `virtual vtkImageData* vtkImageColorViewer::GetInput ( )` [virtual]

```

27.358.4.6  virtual int vtkImageColorViewer::GetOffScreenRendering ( ) [virtual]

27.358.4.7  double vtkImageColorViewer::GetOverlayVisibility ( )

27.358.4.8  virtual int* vtkImageColorViewer::GetPosition ( ) [virtual]

27.358.4.9  virtual int* vtkImageColorViewer::GetSize ( ) [virtual]

27.358.4.10 virtual int vtkImageColorViewer::GetSliceMax ( ) [virtual]

27.358.4.11 virtual int vtkImageColorViewer::GetSliceMin ( ) [virtual]

27.358.4.12 virtual void vtkImageColorViewer::GetSliceRange ( int range[2] ) [inline],[virtual]

27.358.4.13 virtual void vtkImageColorViewer::GetSliceRange ( int & min, int & max ) [virtual]

27.358.4.14 virtual int* vtkImageColorViewer::GetSliceRange ( ) [virtual]

27.358.4.15 virtual const char* vtkImageColorViewer::GetWindowName ( ) [virtual]

27.358.4.16 virtual void vtkImageColorViewer::InstallPipeline ( ) [protected],[virtual]

27.358.4.17 static vtkImageColorViewer* vtkImageColorViewer::New ( ) [static]

```

Examples:

[gdcmrtionplan.cxx](#), and [gdcmrtplan.cxx](#).

```

27.358.4.18 void vtkImageColorViewer::PrintSelf ( ostream & os, vtkIndent indent )

27.358.4.19 virtual void vtkImageColorViewer::Render ( void ) [virtual]

```

Examples:

[gdcmrtionplan.cxx](#), and [gdcmrtplan.cxx](#).

```

27.358.4.20 virtual void vtkImageColorViewer::SetColorLevel ( double s ) [virtual]

27.358.4.21 virtual void vtkImageColorViewer::SetColorWindow ( double s ) [virtual]

27.358.4.22 virtual void vtkImageColorViewer::SetDisplayId ( void * a ) [virtual]

27.358.4.23 virtual void vtkImageColorViewer::SetInput ( vtkImageData * in ) [virtual]

```

Examples:

[gdcmrtionplan.cxx](#), and [gdcmrtplan.cxx](#).

27.358.4.24 `virtual void vtkImageColorViewer::SetInputConnection ( vtkAlgorithmOutput * input )` [virtual]

27.358.4.25 `virtual void vtkImageColorViewer::SetOffScreenRendering ( int )` [virtual]

27.358.4.26 `void vtkImageColorViewer::SetOverlayVisibility ( double vis )`

27.358.4.27 `virtual void vtkImageColorViewer::SetParentId ( void * a )` [virtual]

27.358.4.28 `virtual void vtkImageColorViewer::SetPosition ( int a, int b )` [virtual]

27.358.4.29 `virtual void vtkImageColorViewer::SetPosition ( int a[2] )` [inline],[virtual]

References `SetPosition()`.

Referenced by `SetPosition()`.

27.358.4.30 `virtual void vtkImageColorViewer::SetRenderer ( vtkRenderer * arg )` [virtual]

27.358.4.31 `virtual void vtkImageColorViewer::SetRenderWindow ( vtkRenderWindow * arg )` [virtual]

27.358.4.32 `virtual void vtkImageColorViewer::SetSize ( int a, int b )` [virtual]

Examples:

[gdcmrtionplan.cxx](#), and [gdcmrtplan.cxx](#).

27.358.4.33 `virtual void vtkImageColorViewer::SetSize ( int a[2] )` [inline],[virtual]

References `SetSize()`.

Referenced by `SetSize()`.

27.358.4.34 `virtual void vtkImageColorViewer::SetSlice ( int s )` [virtual]

27.358.4.35 `virtual void vtkImageColorViewer::SetSliceOrientation ( int orientation )` [virtual]

27.358.4.36 `virtual void vtkImageColorViewer::SetSliceOrientationToXY ( )` [inline],[virtual]

References `SLICE_ORIENTATION_XY`.

27.358.4.37 `virtual void vtkImageColorViewer::SetSliceOrientationToXZ ( )` [inline],[virtual]

References `SLICE_ORIENTATION_XZ`.

27.358.4.38 `virtual void vtkImageColorViewer::SetSliceOrientationToYZ ( )` [inline],[virtual]

References `SLICE_ORIENTATION_YZ`.

27.358.4.39 `virtual void vtkImageColorViewer::SetupInteractor ( vtkRenderWindowInteractor * )` [virtual]

Examples:

[gdcmrtonplan.cxx](#), and [gdcmrtpplan.cxx](#).

27.358.4.40 `virtual void vtkImageColorViewer::SetWindowId ( void * a )` [virtual]

27.358.4.41 `virtual void vtkImageColorViewer::UnInstallPipeline ( )` [protected],[virtual]

27.358.4.42 `virtual void vtkImageColorViewer::UpdateDisplayExtent ( )` [virtual]

27.358.4.43 `virtual void vtkImageColorViewer::UpdateOrientation ( )` [protected],[virtual]

27.358.4.44 `vtkImageColorViewer::VTK_LEGACY ( int GetWholeZMin() )`

27.358.4.45 `vtkImageColorViewer::VTK_LEGACY ( int GetWholeZMax() )`

27.358.4.46 `vtkImageColorViewer::VTK_LEGACY ( int GetZSlice() )`

27.358.4.47 `vtkImageColorViewer::VTK_LEGACY ( void SetZSliceint )`

27.358.4.48 `vtkImageColorViewer::vtkBooleanMacro ( OffScreenRendering , int )`

27.358.4.49 `vtkImageColorViewer::vtkGetMacro ( SliceOrientation , int )`

27.358.4.50 `vtkImageColorViewer::vtkGetMacro ( Slice , int )`

27.358.4.51 `vtkImageColorViewer::vtkGetObjectMacro ( RenderWindow , vtkRenderWindow )`

27.358.4.52 `vtkImageColorViewer::vtkGetObjectMacro ( Renderer , vtkRenderer )`

27.358.4.53 `vtkImageColorViewer::vtkGetObjectMacro ( ImageActor , vtkImageActor )`

27.358.4.54 `vtkImageColorViewer::vtkGetObjectMacro ( WindowLevel , vtkImageMapToWindowLevelColors2 )`

27.358.4.55 `vtkImageColorViewer::vtkGetObjectMacro ( InteractorStyle , vtkInteractorStyleImage )`

27.358.4.56 `vtkImageColorViewer::vtkTypeRevisionMacro ( vtkImageColorViewer , vtkObject )`

## 27.358.5 Friends And Related Function Documentation

27.358.5.1 `friend class vtkImageColorViewerCallback` [friend]

## 27.358.6 Member Data Documentation

27.358.6.1 `int vtkImageColorViewer::FirstRender` [protected]

27.358.6.2 `vtkImageActor* vtkImageColorViewer::ImageActor` [protected]

27.358.6.3 `vtkRenderWindowInteractor* vtkImageColorViewer::Interactor` [protected]

27.358.6.4 `vtkInteractorStyleImage*` `vtkImageColorViewer::InteractorStyle` `[protected]`

27.358.6.5 `vtkImageActor*` `vtkImageColorViewer::OverlayImageActor` `[protected]`

27.358.6.6 `vtkRenderer*` `vtkImageColorViewer::Renderer` `[protected]`

27.358.6.7 `vtkRenderWindow*` `vtkImageColorViewer::RenderWindow` `[protected]`

27.358.6.8 `int` `vtkImageColorViewer::Slice` `[protected]`

27.358.6.9 `int` `vtkImageColorViewer::SliceOrientation` `[protected]`

27.358.6.10 `vtkImageMapToWindowLevelColors2*` `vtkImageColorViewer::WindowLevel` `[protected]`


The documentation for this class was generated from the following file:

- [vtkImageColorViewer.h](#)


## 27.359 vtkImageMapToColors16 Class Reference

```
#include <vtkImageMapToColors16.h>
```

Inheritance diagram for `vtkImageMapToColors16`:

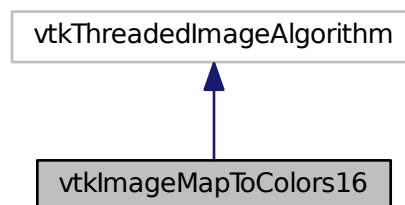

Collaboration diagram for vtkImageMapToColors16:

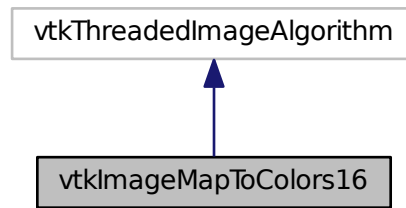

### Public Member Functions

- virtual unsigned long [GetMTime](#) ()
- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- virtual void [SetLookupTable](#) (vtkScalarsToColors \*)
- void [SetOutputFormatToLuminance](#) ()
- void [SetOutputFormatToLuminanceAlpha](#) ()
- void [SetOutputFormatToRGB](#) ()
- void [SetOutputFormatToRGBA](#) ()
- [vtkBooleanMacro](#) (PassAlphaToOutput, int)
- [vtkGetMacro](#) (OutputFormat, int)
- [vtkGetMacro](#) (ActiveComponent, int)
- [vtkGetMacro](#) (PassAlphaToOutput, int)
- [vtkGetObjectMacro](#) (LookupTable, vtkScalarsToColors)
- [vtkSetMacro](#) (OutputFormat, int)
- [vtkSetMacro](#) (ActiveComponent, int)
- [vtkSetMacro](#) (PassAlphaToOutput, int)
- [vtkTypeRevisionMacro](#) (vtkImageMapToColors16, vtkThreadedImageAlgorithm)

### Static Public Member Functions

- static [vtkImageMapToColors16 \\* New](#) ()

### Protected Member Functions

- [vtkImageMapToColors16](#) ()
- [~vtkImageMapToColors16](#) ()
- virtual int [RequestData](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector)
- virtual int [RequestInformation](#) (vtkInformation \*, vtkInformationVector \*\*, vtkInformationVector \*)
- void [ThreadedRequestData](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector, vtkImageData \*\*\*inData, vtkImageData \*\*outData, int extent[6], int id)



## Protected Attributes

- int [ActiveComponent](#)
- int [DataWasPassed](#)
- vtkScalarsToColors \* [LookupTable](#)
- int [OutputFormat](#)
- int [PassAlphaToOutput](#)

## 27.359.1 Constructor & Destructor Documentation

27.359.1.1 `vtkImageMapToColors16::vtkImageMapToColors16 ( )` `[protected]`

27.359.1.2 `vtkImageMapToColors16::~~vtkImageMapToColors16 ( )` `[protected]`

## 27.359.2 Member Function Documentation

27.359.2.1 `virtual unsigned long vtkImageMapToColors16::GetMTime ( )` `[virtual]`

27.359.2.2 `static vtkImageMapToColors16* vtkImageMapToColors16::New ( )` `[static]`

27.359.2.3 `void vtkImageMapToColors16::PrintSelf ( ostream & os, vtkIndent indent )`

27.359.2.4 `virtual int vtkImageMapToColors16::RequestData ( vtkInformation * request, vtkInformationVector ** inputVector, vtkInformationVector * outputVector )` `[protected]`, `[virtual]`

27.359.2.5 `virtual int vtkImageMapToColors16::RequestInformation ( vtkInformation * , vtkInformationVector ** , vtkInformationVector * )` `[protected]`, `[virtual]`

27.359.2.6 `virtual void vtkImageMapToColors16::SetLookupTable ( vtkScalarsToColors * )` `[virtual]`

27.359.2.7 `void vtkImageMapToColors16::SetOutputFormatToLuminance ( )` `[inline]`

27.359.2.8 `void vtkImageMapToColors16::SetOutputFormatToLuminanceAlpha ( )` `[inline]`

27.359.2.9 `void vtkImageMapToColors16::SetOutputFormatToRGB ( )` `[inline]`

27.359.2.10 `void vtkImageMapToColors16::SetOutputFormatToRGBA ( )` `[inline]`

27.359.2.11 `void vtkImageMapToColors16::ThreadedRequestData ( vtkInformation * request, vtkInformationVector ** inputVector, vtkInformationVector * outputVector, vtkImageData *** inData, vtkImageData ** outData, int extent[6], int id )` `[protected]`

27.359.2.12 `vtkImageMapToColors16::vtkBooleanMacro ( PassAlphaToOutput , int )`

27.359.2.13 `vtkImageMapToColors16::vtkGetMacro ( OutputFormat , int )`

27.359.2.14 `vtkImageMapToColors16::vtkGetMacro ( ActiveComponent , int )`

27.359.2.15 `vtkImageMapToColors16::vtkGetMacro ( PassAlphaToOutput , int )`

27.359.2.16 `vtkImageMapToColors16::vtkGetObjectMacro ( LookupTable , vtkScalarsToColors )`

27.359.2.17 `vtkImageMapToColors16::vtkSetMacro ( OutputFormat , int )`

27.359.2.18 `vtkImageMapToColors16::vtkSetMacro ( ActiveComponent , int )`

27.359.2.19 `vtkImageMapToColors16::vtkSetMacro ( PassAlphaToOutput , int )`

27.359.2.20 `vtkImageMapToColors16::vtkTypeRevisionMacro ( vtkImageMapToColors16 , vtkThreadedImageAlgorithm )`

### 27.359.3 Member Data Documentation

27.359.3.1 `int vtkImageMapToColors16::ActiveComponent` [protected]

27.359.3.2 `int vtkImageMapToColors16::DataWasPassed` [protected]

27.359.3.3 `vtkScalarsToColors* vtkImageMapToColors16::LookupTable` [protected]

27.359.3.4 `int vtkImageMapToColors16::OutputFormat` [protected]

27.359.3.5 `int vtkImageMapToColors16::PassAlphaToOutput` [protected]


The documentation for this class was generated from the following file:

- [vtkImageMapToColors16.h](#)


## 27.360 vtkImageMapToWindowLevelColors2 Class Reference

```
#include <vtkImageMapToWindowLevelColors2.h>
```

Inheritance diagram for `vtkImageMapToWindowLevelColors2`:

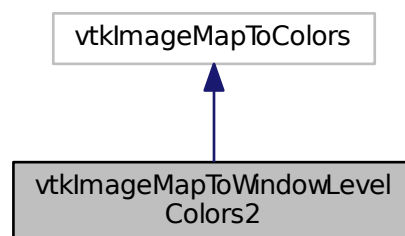

Collaboration diagram for vtkImageMapToWindowLevelColors2:

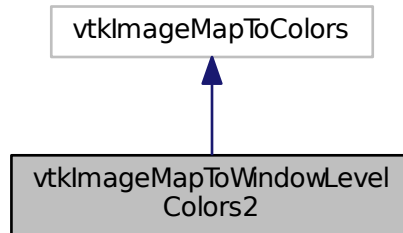

### Public Member Functions

- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkGetMacro](#) (Window, double)
- [vtkGetMacro](#) (Level, double)
- [vtkSetMacro](#) (Window, double)
- [vtkSetMacro](#) (Level, double)
- [vtkTypeRevisionMacro](#) (vtkImageMapToWindowLevelColors2, vtkImageMapToColors)

### Static Public Member Functions

- static [vtkImageMapToWindowLevelColors2 \\* New](#) ()

### Protected Member Functions

- [vtkImageMapToWindowLevelColors2](#) ()
- [~vtkImageMapToWindowLevelColors2](#) ()
- virtual int [RequestData](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector)
- virtual int [RequestInformation](#) (vtkInformation \*, vtkInformationVector \*\*, vtkInformationVector \*)
- void [ThreadedRequestData](#) (vtkInformation \*request, vtkInformationVector \*\*inputVector, vtkInformationVector \*outputVector, vtkImageData \*\*\*inData, vtkImageData \*\*outData, int extent[6], int id)

### Protected Attributes

- double [Level](#)
- double [Window](#)

### 27.360.1 Constructor & Destructor Documentation

27.360.1.1 `vtkImageMapToWindowLevelColors2::vtkImageMapToWindowLevelColors2 ( )` [protected]

27.360.1.2 `vtkImageMapToWindowLevelColors2::~~vtkImageMapToWindowLevelColors2 ( )` [protected]

### 27.360.2 Member Function Documentation

27.360.2.1 `static vtkImageMapToWindowLevelColors2* vtkImageMapToWindowLevelColors2::New ( )` [static]

27.360.2.2 `void vtkImageMapToWindowLevelColors2::PrintSelf ( ostream & os, vtkIndent indent )`

27.360.2.3 `virtual int vtkImageMapToWindowLevelColors2::RequestData ( vtkInformation * request, vtkInformationVector ** inputVector, vtkInformationVector * outputVector )` [protected],[virtual]

27.360.2.4 `virtual int vtkImageMapToWindowLevelColors2::RequestInformation ( vtkInformation * , vtkInformationVector ** , vtkInformationVector * )` [protected],[virtual]

27.360.2.5 `void vtkImageMapToWindowLevelColors2::ThreadedRequestData ( vtkInformation * request, vtkInformationVector ** inputVector, vtkInformationVector * outputVector, vtkImageData *** inData, vtkImageData ** outData, int extent[6], int id )` [protected]

27.360.2.6 `vtkImageMapToWindowLevelColors2::vtkGetMacro ( Window , double )`

27.360.2.7 `vtkImageMapToWindowLevelColors2::vtkGetMacro ( Level , double )`

27.360.2.8 `vtkImageMapToWindowLevelColors2::vtkSetMacro ( Window , double )`

27.360.2.9 `vtkImageMapToWindowLevelColors2::vtkSetMacro ( Level , double )`

27.360.2.10 `vtkImageMapToWindowLevelColors2::vtkTypeRevisionMacro ( vtkImageMapToWindowLevelColors2 , vtkImageMapToColors )`

### 27.360.3 Member Data Documentation

27.360.3.1 `double vtkImageMapToWindowLevelColors2::Level` [protected]

27.360.3.2 `double vtkImageMapToWindowLevelColors2::Window` [protected]


The documentation for this class was generated from the following file:

- [vtkImageMapToWindowLevelColors2.h](#)


## 27.361 vtkImagePlanarComponentsToComponents Class Reference

```
#include <vtkImagePlanarComponentsToComponents.h>
```

Inheritance diagram for vtkImagePlanarComponentsToComponents:

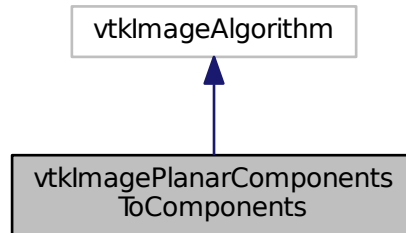

Collaboration diagram for vtkImagePlanarComponentsToComponents:

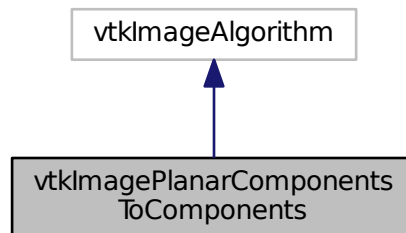

### Public Member Functions

- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkTypeRevisionMacro](#) ([vtkImagePlanarComponentsToComponents](#), vtkImageAlgorithm)

### Static Public Member Functions

- static [vtkImagePlanarComponentsToComponents](#) \* [New](#) ()

### Protected Member Functions

- [vtkImagePlanarComponentsToComponents](#) ()
- [~vtkImagePlanarComponentsToComponents](#) ()
- virtual int [RequestData](#) (vtkInformation \*, vtkInformationVector \*\*, vtkInformationVector \*)

### 27.361.1 Constructor & Destructor Documentation

27.361.1.1 `vtkImagePlanarComponentsToComponents::vtkImagePlanarComponentsToComponents ( )` [protected]

27.361.1.2 `vtkImagePlanarComponentsToComponents::~~vtkImagePlanarComponentsToComponents ( )` [inline], [protected]

### 27.361.2 Member Function Documentation

27.361.2.1 `static vtkImagePlanarComponentsToComponents* vtkImagePlanarComponentsToComponents::New ( )` [static]

27.361.2.2 `void vtkImagePlanarComponentsToComponents::PrintSelf ( ostream & os, vtkIndent indent )`

27.361.2.3 `virtual int vtkImagePlanarComponentsToComponents::RequestData ( vtkInformation *, vtkInformationVector **, vtkInformationVector * )` [protected], [virtual]

27.361.2.4 `vtkImagePlanarComponentsToComponents::vtkTypeRevisionMacro ( vtkImagePlanarComponentsToComponents, vtkImageAlgorithm )`


The documentation for this class was generated from the following file:

- [vtkImagePlanarComponentsToComponents.h](#)


## 27.362 vtkImageRGBToYBR Class Reference

```
#include <vtkImageRGBToYBR.h>
```

Inheritance diagram for `vtkImageRGBToYBR`:

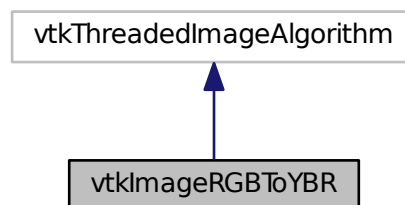

Collaboration diagram for vtkImageRGBToYBR:

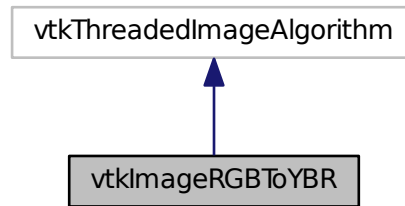

### Public Member Functions

- void `PrintSelf` (ostream &os, vtkIndent indent)
- `vtkTypeRevisionMacro` (`vtkImageRGBToYBR`, `vtkThreadedImageAlgorithm`)

### Static Public Member Functions

- static `vtkImageRGBToYBR * New` ()

### Protected Member Functions

- `vtkImageRGBToYBR` ()
- `~vtkImageRGBToYBR` ()
- void `ThreadedExecute` (vtkImageData \*inData, vtkImageData \*outData, int ext[6], int id)

## 27.362.1 Constructor & Destructor Documentation

27.362.1.1 `vtkImageRGBToYBR::vtkImageRGBToYBR ( )` [protected]

27.362.1.2 `vtkImageRGBToYBR::~~vtkImageRGBToYBR ( )` [inline], [protected]

## 27.362.2 Member Function Documentation

27.362.2.1 `static vtkImageRGBToYBR* vtkImageRGBToYBR::New ( )` [static]

27.362.2.2 `void vtkImageRGBToYBR::PrintSelf ( ostream & os, vtkIndent indent )`

27.362.2.3 `void vtkImageRGBToYBR::ThreadedExecute ( vtkImageData * inData, vtkImageData * outData, int ext[6], int id )`  
[protected]

27.362.2.4 `vtkImageRGBToYBR::vtkTypeRevisionMacro ( vtkImageRGBToYBR , vtkThreadedImageAlgorithm )`


The documentation for this class was generated from the following file:

- [vtkImageRGBToYBR.h](#)


## 27.363 vtkImageYBRToRGB Class Reference

```
#include <vtkImageYBRToRGB.h>
```

Inheritance diagram for vtkImageYBRToRGB:

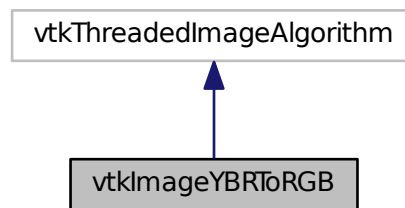

Collaboration diagram for vtkImageYBRToRGB:

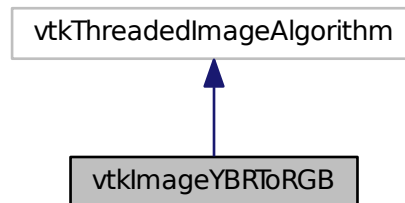

### Public Member Functions

- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkTypeRevisionMacro](#) (vtkImageYBRToRGB, vtkThreadedImageAlgorithm)

### Static Public Member Functions

- static [vtkImageYBRToRGB \\* New](#) ()



## Protected Member Functions

- [vtkImageYBRToRGB](#) ()
- [~vtkImageYBRToRGB](#) ()
- void [ThreadedExecute](#) (vtkImageData \*inData, vtkImageData \*outData, int ext[6], int id)

## 27.363.1 Constructor & Destructor Documentation

27.363.1.1 `vtkImageYBRToRGB::vtkImageYBRToRGB ( )` [protected]

27.363.1.2 `vtkImageYBRToRGB::~~vtkImageYBRToRGB ( )` [inline], [protected]

## 27.363.2 Member Function Documentation

27.363.2.1 `static vtkImageYBRToRGB* vtkImageYBRToRGB::New ( )` [static]

27.363.2.2 `void vtkImageYBRToRGB::PrintSelf ( ostream & os, vtkIndent indent )`

27.363.2.3 `void vtkImageYBRToRGB::ThreadedExecute ( vtkImageData * inData, vtkImageData * outData, int ext[6], int id )`  
[protected]

27.363.2.4 `vtkImageYBRToRGB::vtkTypeRevisionMacro ( vtkImageYBRToRGB , vtkThreadedImageAlgorithm )`


The documentation for this class was generated from the following file:

- [vtkImageYBRToRGB.h](#)


## 27.364 vtkLookupTable16 Class Reference

```
#include <vtkLookupTable16.h>
```

Inheritance diagram for vtkLookupTable16:

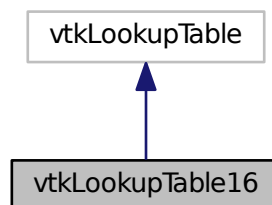

Collaboration diagram for vtkLookupTable16:

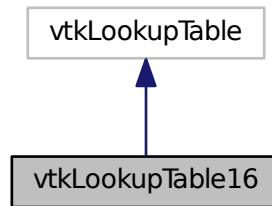

### Public Member Functions

- void [Build](#) ()
- unsigned short \* [GetPointer](#) (const vtkIdType id)
- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- void [SetNumberOfTableValues](#) (vtkIdType number)
- [vtkTypeRevisionMacro](#) ([vtkLookupTable16](#), vtkLookupTable)
- unsigned char \* [WritePointer](#) (const vtkIdType id, const int number)

### Static Public Member Functions

- static [vtkLookupTable16](#) \* [New](#) ()

### Protected Member Functions

- [vtkLookupTable16](#) (int size=256, int ext=256)
- [~vtkLookupTable16](#) ()
- void [MapScalarsThroughTable2](#) (void \*input, unsigned char \*output, int inputDataType, int numberOfValues, int inputIncrement, int outputFormat)

### Protected Attributes

- vtkUnsignedShortArray \* [Table16](#)

## 27.364.1 Constructor & Destructor Documentation

27.364.1.1 `vtkLookupTable16::vtkLookupTable16 ( int size = 256, int ext = 256 )` [protected]

27.364.1.2 `vtkLookupTable16::~~vtkLookupTable16 ( )` [protected]

## 27.364.2 Member Function Documentation

- 27.364.2.1 void vtkLookupTable16::Build ( )
- 27.364.2.2 unsigned short\* vtkLookupTable16::GetPointer ( const vtkIdType *id* ) [inline]
- 27.364.2.3 void vtkLookupTable16::MapScalarsThroughTable2 ( void \* *input*, unsigned char \* *output*, int *inputDataType*, int *numberOfValues*, int *inputIncrement*, int *outputFormat* ) [protected]
- 27.364.2.4 static vtkLookupTable16\* vtkLookupTable16::New ( ) [static]
- 27.364.2.5 void vtkLookupTable16::PrintSelf ( ostream & *os*, vtkIndent *indent* )
- 27.364.2.6 void vtkLookupTable16::SetNumberOfTableValues ( vtkIdType *number* )
- 27.364.2.7 vtkLookupTable16::vtkTypeRevisionMacro ( vtkLookupTable16 , vtkLookupTable )
- 27.364.2.8 unsigned char \* vtkLookupTable16::WritePointer ( const vtkIdType *id*, const int *number* ) [inline]

References Table16.

### 27.364.3 Member Data Documentation

- 27.364.3.1 vtkUnsignedShortArray\* vtkLookupTable16::Table16 [protected]

Referenced by WritePointer().


The documentation for this class was generated from the following file:

- [vtkLookupTable16.h](#)


## 27.365 vtkRTStructSetProperties Class Reference

```
#include <vtkRTStructSetProperties.h>
```

Inheritance diagram for vtkRTStructSetProperties:

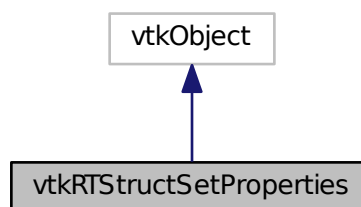

Collaboration diagram for vtkRTStructSetProperties:

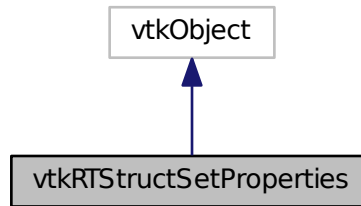

## Public Member Functions

- void [AddContourReferencedFrameOfReference](#) (vtkIdType pdnum, const char \*classuid, const char \*instanceuid)
- void [AddReferencedFrameOfReference](#) (const char \*classuid, const char \*instanceuid)
- void [AddStructureSetROI](#) (int roinumber, const char \*reframerefid, const char \*roiname, const char \*ROI←  
GenerationAlgorithm, const char \*ROIDescription=0)
- void [AddStructureSetROIObservation](#) (int refnumber, int observationnumber, const char \*rtroiinterpretedtype,  
const char \*roiinterpreter, const char \*roiobservationlabel=0)
- virtual void [Clear](#) ()
- virtual void [DeepCopy](#) (vtkRTStructSetProperties \*p)
- const char \* [GetContourReferencedFrameOfReferenceClassUID](#) (vtkIdType pdnum, vtkIdType id)
- const char \* [GetContourReferencedFrameOfReferenceInstanceUID](#) (vtkIdType pdnum, vtkIdType id)
- vtkIdType [GetNumberOfContourReferencedFrameOfReferences](#) ()
- vtkIdType [GetNumberOfContourReferencedFrameOfReferences](#) (vtkIdType pdnum)
- vtkIdType [GetNumberOfReferencedFrameOfReferences](#) ()
- vtkIdType [GetNumberOfStructureSetROIs](#) ()
- const char \* [GetReferencedFrameOfReferenceClassUID](#) (vtkIdType id)
- const char \* [GetReferencedFrameOfReferenceInstanceUID](#) (vtkIdType id)
- int [GetStructureSetObservationNumber](#) (vtkIdType id)
- const char \* [GetStructureSetROIDescription](#) (vtkIdType id)
- const char \* [GetStructureSetROIGenerationAlgorithm](#) (vtkIdType)
- const char \* [GetStructureSetROIName](#) (vtkIdType)
- int [GetStructureSetROINumber](#) (vtkIdType id)
- const char \* [GetStructureSetROIObservationLabel](#) (vtkIdType id)
- const char \* [GetStructureSetROIRefFrameRefUID](#) (vtkIdType)
- const char \* [GetStructureSetRTROIInterpretedType](#) (vtkIdType id)
- void [PrintSelf](#) (ostream &os, vtkIndent indent)
- [vtkGetStringMacro](#) (StructureSetLabel)
- [vtkGetStringMacro](#) (StructureSetName)
- [vtkGetStringMacro](#) (StructureSetDate)
- [vtkGetStringMacro](#) (StructureSetTime)
- [vtkGetStringMacro](#) (SOPInstanceUID)
- [vtkGetStringMacro](#) (StudyInstanceUID)
- [vtkGetStringMacro](#) (SeriesInstanceUID)

- [vtkGetStringMacro \(ReferenceSeriesInstanceUID\)](#)
- [vtkGetStringMacro \(ReferenceFrameOfReferenceUID\)](#)
- [vtkSetStringMacro \(StructureSetLabel\)](#)
- [vtkSetStringMacro \(StructureSetName\)](#)
- [vtkSetStringMacro \(StructureSetDate\)](#)
- [vtkSetStringMacro \(StructureSetTime\)](#)
- [vtkSetStringMacro \(SOPInstanceUID\)](#)
- [vtkSetStringMacro \(StudyInstanceUID\)](#)
- [vtkSetStringMacro \(SeriesInstanceUID\)](#)
- [vtkSetStringMacro \(ReferenceSeriesInstanceUID\)](#)
- [vtkSetStringMacro \(ReferenceFrameOfReferenceUID\)](#)
- [vtkTypeRevisionMacro \(vtkRTStructSetProperties, vtkObject\)](#)

### Static Public Member Functions

- static [vtkRTStructSetProperties \\* New \(\)](#)

### Protected Member Functions

- [vtkRTStructSetProperties \(\)](#)
- [~vtkRTStructSetProperties \(\)](#)

### Protected Attributes

- [vtkRTStructSetPropertiesInternals \\* Internals](#)
- [char \\* ReferenceFrameOfReferenceUID](#)
- [char \\* ReferenceSeriesInstanceUID](#)
- [char \\* SeriesInstanceUID](#)
- [char \\* SOPInstanceUID](#)
- [char \\* StructureSetDate](#)
- [char \\* StructureSetLabel](#)
- [char \\* StructureSetName](#)
- [char \\* StructureSetTime](#)
- [char \\* StudyInstanceUID](#)

## 27.365.1 Detailed Description

Examples:

[GenerateRTSTRUCT.cxx](#).

## 27.365.2 Constructor & Destructor Documentation

27.365.2.1 [vtkRTStructSetProperties::vtkRTStructSetProperties \( \)](#) [protected]

27.365.2.2 [vtkRTStructSetProperties::~~vtkRTStructSetProperties \( \)](#) [protected]

## 27.365.3 Member Function Documentation

- 27.365.3.1 void vtkRTStructSetProperties::AddContourReferencedFrameOfReference ( vtkIdType *pdnum*, const char \* *classuid*, const char \* *instanceuid* )
- 27.365.3.2 void vtkRTStructSetProperties::AddReferencedFrameOfReference ( const char \* *classuid*, const char \* *instanceuid* )
- 27.365.3.3 void vtkRTStructSetProperties::AddStructureSetROI ( int *roinumber*, const char \* *refframerefid*, const char \* *roiname*, const char \* *ROIGenerationAlgorithm*, const char \* *ROIDescription* = 0 )
- 27.365.3.4 void vtkRTStructSetProperties::AddStructureSetROIObservation ( int *refnumber*, int *observationnumber*, const char \* *rtroiinterpretedtype*, const char \* *roiinterpreter*, const char \* *roiobservationlabel* = 0 )
- 27.365.3.5 virtual void vtkRTStructSetProperties::Clear ( ) [virtual]
- 27.365.3.6 virtual void vtkRTStructSetProperties::DeepCopy ( vtkRTStructSetProperties \* *p* ) [virtual]
- 27.365.3.7 const char\* vtkRTStructSetProperties::GetContourReferencedFrameOfReferenceClassUID ( vtkIdType *pdnum*, vtkIdType *id* )
- 27.365.3.8 const char\* vtkRTStructSetProperties::GetContourReferencedFrameOfReferenceInstanceUID ( vtkIdType *pdnum*, vtkIdType *id* )
- 27.365.3.9 vtkIdType vtkRTStructSetProperties::GetNumberOfContourReferencedFrameOfReferences ( )
- 27.365.3.10 vtkIdType vtkRTStructSetProperties::GetNumberOfContourReferencedFrameOfReferences ( vtkIdType *pdnum* )
- 27.365.3.11 vtkIdType vtkRTStructSetProperties::GetNumberOfReferencedFrameOfReferences ( )
- 27.365.3.12 vtkIdType vtkRTStructSetProperties::GetNumberOfStructureSetROIs ( )
- 27.365.3.13 const char\* vtkRTStructSetProperties::GetReferencedFrameOfReferenceClassUID ( vtkIdType *id* )
- 27.365.3.14 const char\* vtkRTStructSetProperties::GetReferencedFrameOfReferenceInstanceUID ( vtkIdType *id* )
- 27.365.3.15 int vtkRTStructSetProperties::GetStructureSetObservationNumber ( vtkIdType *id* )
- 27.365.3.16 const char\* vtkRTStructSetProperties::GetStructureSetROIDescription ( vtkIdType *id* )
- 27.365.3.17 const char\* vtkRTStructSetProperties::GetStructureSetROIGenerationAlgorithm ( vtkIdType )
- 27.365.3.18 const char\* vtkRTStructSetProperties::GetStructureSetROIName ( vtkIdType )
- 27.365.3.19 int vtkRTStructSetProperties::GetStructureSetROINumber ( vtkIdType *id* )
- 27.365.3.20 const char\* vtkRTStructSetProperties::GetStructureSetROIObservationLabel ( vtkIdType *id* )
- 27.365.3.21 const char\* vtkRTStructSetProperties::GetStructureSetROIRefFrameRefUID ( vtkIdType )
- 27.365.3.22 const char\* vtkRTStructSetProperties::GetStructureSetRTROIInterpretedType ( vtkIdType *id* )

27.365.3.23 static vtkRTStructSetProperties\* vtkRTStructSetProperties::New ( ) [static]

Examples:

[GenerateRTSTRUCT.cxx](#).

27.365.3.24 void vtkRTStructSetProperties::PrintSelf ( ostream & os, vtkIndent indent )

27.365.3.25 vtkRTStructSetProperties::vtkGetStringMacro ( StructureSetLabel )

27.365.3.26 vtkRTStructSetProperties::vtkGetStringMacro ( StructureSetName )

27.365.3.27 vtkRTStructSetProperties::vtkGetStringMacro ( StructureSetDate )

27.365.3.28 vtkRTStructSetProperties::vtkGetStringMacro ( StructureSetTime )

27.365.3.29 vtkRTStructSetProperties::vtkGetStringMacro ( SOPInstanceUID )

27.365.3.30 vtkRTStructSetProperties::vtkGetStringMacro ( StudyInstanceUID )

27.365.3.31 vtkRTStructSetProperties::vtkGetStringMacro ( SeriesInstanceUID )

27.365.3.32 vtkRTStructSetProperties::vtkGetStringMacro ( ReferenceSeriesInstanceUID )

27.365.3.33 vtkRTStructSetProperties::vtkGetStringMacro ( ReferenceFrameOfReferenceUID )

27.365.3.34 vtkRTStructSetProperties::vtkSetStringMacro ( StructureSetLabel )

27.365.3.35 vtkRTStructSetProperties::vtkSetStringMacro ( StructureSetName )

27.365.3.36 vtkRTStructSetProperties::vtkSetStringMacro ( StructureSetDate )

27.365.3.37 vtkRTStructSetProperties::vtkSetStringMacro ( StructureSetTime )

27.365.3.38 vtkRTStructSetProperties::vtkSetStringMacro ( SOPInstanceUID )

27.365.3.39 vtkRTStructSetProperties::vtkSetStringMacro ( StudyInstanceUID )

27.365.3.40 vtkRTStructSetProperties::vtkSetStringMacro ( SeriesInstanceUID )

27.365.3.41 vtkRTStructSetProperties::vtkSetStringMacro ( ReferenceSeriesInstanceUID )

27.365.3.42 vtkRTStructSetProperties::vtkSetStringMacro ( ReferenceFrameOfReferenceUID )

27.365.3.43 vtkRTStructSetProperties::vtkTypeRevisionMacro ( vtkRTStructSetProperties , vtkObject )

## 27.365.4 Member Data Documentation

27.365.4.1 vtkRTStructSetPropertiesInternals\* vtkRTStructSetProperties::Internals [protected]

27.365.4.2 char\* vtkRTStructSetProperties::ReferenceFrameOfReferenceUID [protected]

27.365.4.3 `char* vtkRTStructSetProperties::ReferenceSeriesInstanceUID` [protected]

27.365.4.4 `char* vtkRTStructSetProperties::SeriesInstanceUID` [protected]

27.365.4.5 `char* vtkRTStructSetProperties::SOPInstanceUID` [protected]

27.365.4.6 `char* vtkRTStructSetProperties::StructureSetDate` [protected]

27.365.4.7 `char* vtkRTStructSetProperties::StructureSetLabel` [protected]

27.365.4.8 `char* vtkRTStructSetProperties::StructureSetName` [protected]

27.365.4.9 `char* vtkRTStructSetProperties::StructureSetTime` [protected]

27.365.4.10 `char* vtkRTStructSetProperties::StudyInstanceUID` [protected]

The documentation for this class was generated from the following file:

- [vtkRTStructSetProperties.h](#)

## 27.366 gdcm::Waveform Class Reference

[Waveform](#) class.

```
#include <gdcmWaveform.h>
```

### Public Member Functions

- [Waveform](#) ()

### 27.366.1 Detailed Description

[Waveform](#) class.

### 27.366.2 Constructor & Destructor Documentation

27.366.2.1 `gdcm::Waveform::Waveform ( )` [inline]

The documentation for this class was generated from the following file:

- [gdcmWaveform.h](#)


## 27.367 gdcm::Writer Class Reference

[Writer](#) ala DOM (Document [Object](#) Model) This class is a non-validating writer, it will only performs well- formedness check only.


```
#include <gdcmWriter.h>
```



Inheritance diagram for gdcm::Writer:

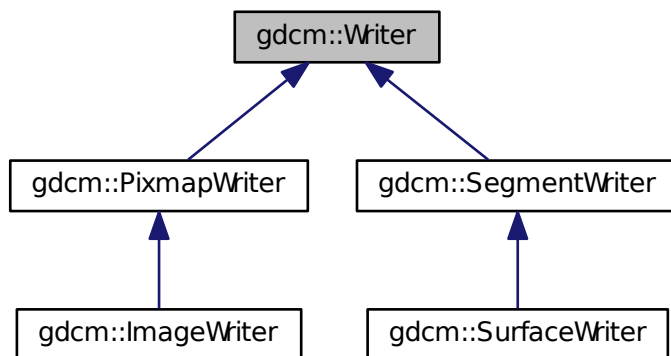

Collaboration diagram for `gdcm::Writer`:

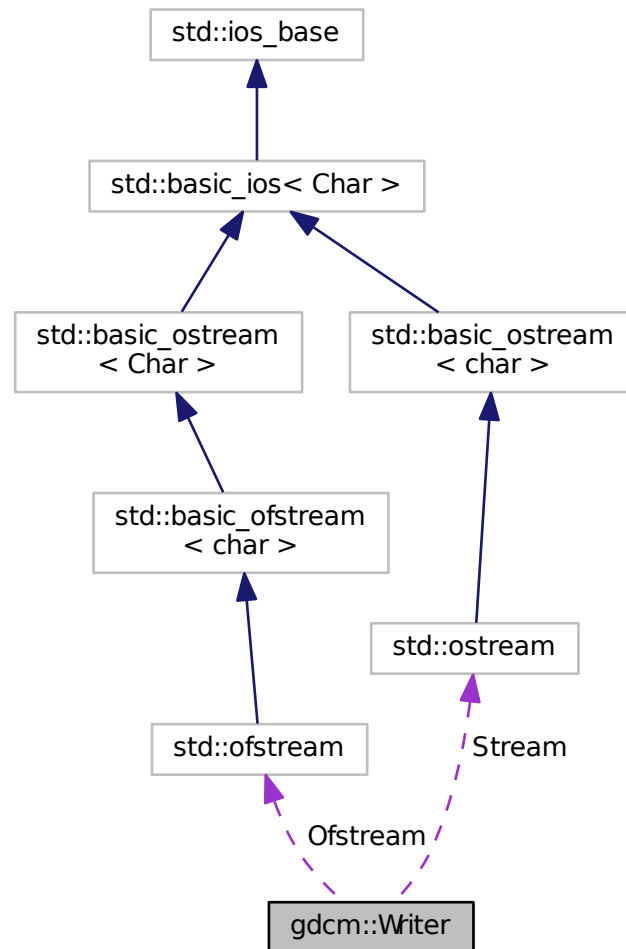

## Public Member Functions

- [Writer](#) ()
- virtual [~Writer](#) ()
- void [CheckFileMetaInformationOff](#) ()
- void [CheckFileMetaInformationOn](#) ()
- [File](#) & [GetFile](#) ()
- void [SetCheckFileMetaInformation](#) (bool b)  
*Undocumented function, do not use (= leave default)*
- void [SetFile](#) (const [File](#) &f)  
*Set/Get the DICOM file ([DataSet](#) + Header)*
- void [SetFileName](#) (const char \*filename\_native)

*Set the filename of DICOM file to write:*

- void [SetStream](#) (std::ostream &output\_stream)

*Set user ostream buffer.*

- virtual bool [Write](#) ()

*Main function to tell the writer to write.*

## Protected Member Functions

- std::ostream \* [GetStreamPtr](#) () const
- void [SetWriteDataSetOnly](#) (bool b)

## Protected Attributes

- std::ofstream \* [Ofstream](#)
- std::ostream \* [Stream](#)

## Friends

- class [StreamImageWriter](#)

## 27.367.1 Detailed Description

[Writer](#) ala DOM (Document [Object](#) Model) This class is a non-validating writer, it will only performs well- formedness check only.

Detailed description here To avoid GDCM being yet another broken DICOM lib we try to be user level and avoid writing illegal stuff (odd length, non-zero value for [Item](#) start/end length ...) Therefore you cannot (well unless you are really smart) write DICOM with even length tag. All the checks are consider basics:

- Correct Meta Information Header (see [gdcm::FileMetaInformation](#))
- Zero value for [Item](#) Length (0xfffe, 0xe00d/0xe0dd)
- Even length for any elements
- Alphabetical order for elements (garanteed by design of internals)
- 32bits [VR](#) will be rewritten with 00

### Warning

[gdcm::Writer](#) cannot write a [DataSet](#) if no SOP Instance UID (0008,0018) is found, unless a [DICOMDIR](#) is being written out

### See also

[Reader DataSet File](#)

### Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [DuplicatePCDE.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenAllVR.cxx](#), [GenFakelIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [HelloWorld.cxx](#), [LargeVRDSExplicit.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [rle2img.cxx](#), and [StreamImageReaderTest.cxx](#).

## 27.367.2 Constructor & Destructor Documentation

27.367.2.1 `gdcm::Writer::Writer ( )`

27.367.2.2 `virtual gdcm::Writer::~~Writer ( )` `[virtual]`

## 27.367.3 Member Function Documentation

27.367.3.1 `void gdcm::Writer::CheckFileMetaInformationOff ( )` `[inline]`

Examples:

[CreateFakeRTDOSE.cxx](#), [FixBrokenJ2K.cxx](#), and [HelloWorld.cxx](#).

27.367.3.2 `void gdcm::Writer::CheckFileMetaInformationOn ( )` `[inline]`

27.367.3.3 `File& gdcm::Writer::GetFile ( )` `[inline]`

Examples:

[CreateJPIPDataSet.cxx](#), [EncapsulateFileInRawData.cxx](#), [Extracting\\_All\\_Resolution.cxx](#), [Fake\\_Image\\_Using\\_Stream\\_Image\\_Writer.cxx](#), [GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), [GetSubSequenceData.cxx](#), [iU22tomultisc.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), [rle2img.cxx](#), and [StreamImageReaderTest.cxx](#).

27.367.3.4 `std::ostream* gdcm::Writer::GetStreamPtr ( ) const` `[inline]`, `[protected]`

27.367.3.5 `void gdcm::Writer::SetCheckFileMetaInformation ( bool b )` `[inline]`

Undocumented function, do not use (= leave default)

Examples:

[GenAllVR.cxx](#), [GenFakeIdentifyFile.cxx](#), and [PatchFile.cxx](#).

27.367.3.6 `void gdcm::Writer::SetFile ( const File & f )` `[inline]`

Set/Get the DICOM file ([DataSet](#) + Header)

Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CompressImage.cxx](#), [CreateFakeRTDOSE.cxx](#), [DuplicatePCDE.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenFakeImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [HelloWorld.cxx](#), [LargeVRDSEExplicit.cxx](#), [MergeTwoFiles.cxx](#), [NewSequence.cs](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), and [rle2img.cxx](#).

27.367.3.7 `void gdcm::Writer::SetFileName ( const char * filename_native )`

Set the filename of DICOM file to write:

## Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CompressImage.cxx](#), [CreateARGBImage.cxx](#), [CreateCMYKImage.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [csa2img.cxx](#), [DuplicatePCDE.cxx](#), [EncapsulateFileInRawData.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenAllVR.cxx](#), [GenFakelIdentifyFile.cxx](#), [GenFakelImage.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [GetSubSequenceData.cxx](#), [HelloVizWorld.cxx](#), [HelloWorld.cxx](#), [iU22tomultisc.cxx](#), [LargeVRDSEExplicit.cxx](#), [MergeTwoFiles.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), and [rle2img.cxx](#).

**27.367.3.8** `void gdcm::Writer::SetStream ( std::ostream & output_stream ) [inline]`

Set user ostream buffer.

**27.367.3.9** `void gdcm::Writer::SetWriteDataSetOnly ( bool b ) [inline], [protected]`

**27.367.3.10** `virtual bool gdcm::Writer::Write ( ) [virtual]`

Main function to tell the writer to write.

Reimplemented in [gdcm::PixmapWriter](#), [gdcm::ImageWriter](#), [gdcm::SurfaceWriter](#), and [gdcm::SegmentWriter](#).

## Examples:

[ChangePrivateTags.cxx](#), [ChangeSequenceUltrasound.cxx](#), [ClinicalTrialAnnotate.cxx](#), [CreateFakeRTDOSE.cxx](#), [CreateJPIPDataSet.cxx](#), [DuplicatePCDE.cxx](#), [EncapsulateFileInRawData.cxx](#), [FixBrokenJ2K.cxx](#), [FixJAIBugJPEGLS.cxx](#), [GenAllVR.cxx](#), [GenFakelIdentifyFile.cxx](#), [GenLongSeqs.cxx](#), [GenSeqs.cxx](#), [HelloWorld.cxx](#), [LargeVRDSEExplicit.cxx](#), [PatchFile.cxx](#), [pmsct\\_rgb1.cxx](#), [QIDO-RS.cxx](#), and [rle2img.cxx](#).

**27.367.4 Friends And Related Function Documentation**

**27.367.4.1** `friend class StreamImageWriter [friend]`

**27.367.5 Member Data Documentation**

**27.367.5.1** `std::ofstream* gdcm::Writer::Ofstream [protected]`

**27.367.5.2** `std::ostream* gdcm::Writer::Stream [protected]`

The documentation for this class was generated from the following file:


- [gdcmWriter.h](#)

**27.368 gdcm::XMLDictReader Class Reference**


Class for representing a [XMLDictReader](#).

```
#include <gdcmXMLDictReader.h>
```

Inheritance diagram for `gdcm::XMLDictReader`:

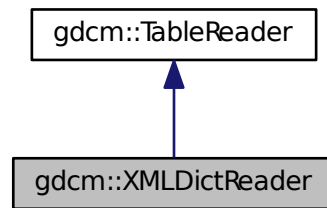

Collaboration diagram for `gdcm::XMLDictReader`:

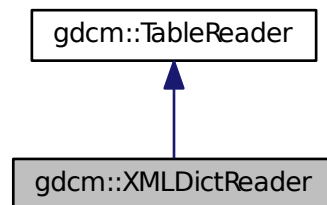

### Public Member Functions

- [XMLDictReader](#) ()
- [~XMLDictReader](#) ()
- void [CharacterDataHandler](#) (const char \*data, int length)
- void [EndElement](#) (const char \*name)
- const [Dict](#) & [GetDict](#) ()
- void [StartElement](#) (const char \*name, const char \*\*atts)

### Protected Member Functions

- void [HandleDescription](#) (const char \*\*atts)
- void [HandleEntry](#) (const char \*\*atts)

## 27.368.1 Detailed Description

Class for representing a [XMLDictReader](#).

## Note

bla Will read the DICOMV3.xml file

## 27.368.2 Constructor & Destructor Documentation

27.368.2.1 `gdcm::XMLDictReader::XMLDictReader ( )`

27.368.2.2 `gdcm::XMLDictReader::~~XMLDictReader ( )` `[inline]`

## 27.368.3 Member Function Documentation

27.368.3.1 `void gdcm::XMLDictReader::CharacterDataHandler ( const char * data, int length )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).

27.368.3.2 `void gdcm::XMLDictReader::EndElement ( const char * name )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).

27.368.3.3 `const Dict& gdcm::XMLDictReader::GetDict ( )` `[inline]`

27.368.3.4 `void gdcm::XMLDictReader::HandleDescription ( const char ** atts )` `[protected]`

27.368.3.5 `void gdcm::XMLDictReader::HandleEntry ( const char ** atts )` `[protected]`

27.368.3.6 `void gdcm::XMLDictReader::StartElement ( const char * name, const char ** atts )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).


The documentation for this class was generated from the following file:

- [gdcmXMLDictReader.h](#)

## 27.369 gdcm::XMLPrinter Class Reference

```
#include <gdcmXMLPrinter.h>
```

Collaboration diagram for `gdcm::XMLPrinter`:

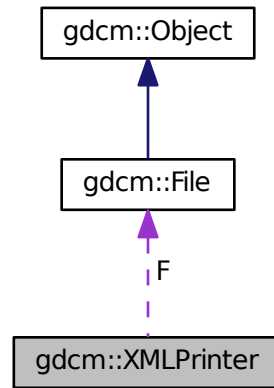

## Public Types

- enum `PrintStyles` {  
`OnlyUUID` = 0,  
`LOADBULKDATA` = 1 }

## Public Member Functions

- `XMLPrinter` ()
- virtual `~XMLPrinter` ()
- `PrintStyles` `GetPrintStyle` () const
- virtual void `HandleBulkData` (const char \*uuid, const `TransferSyntax` &ts, const char \*bulkdata, size\_t bulklen)
- void `Print` (std::ostream &os)
- void `PrintDataSet` (const `DataSet` &ds, const `TransferSyntax` &ts, std::ostream &os)
- void `SetFile` (`File` const &f)
- void `SetStyle` (`PrintStyles` ps)

## Protected Member Functions

- VR `PrintDataElement` (std::ostream &os, const `Dicts` &dicts, const `DataSet` &ds, const `DataElement` &de, const `TransferSyntax` &ts)
- void `PrintSQ` (const `SequenceOfItems` \*sqi, const `TransferSyntax` &ts, std::ostream &os)

## Protected Attributes

- const `File` \* `F`
- `PrintStyles` `PrintStyle`



## 27.369.1 Member Enumeration Documentation

### 27.369.1.1 enum gdcm::XMLPrinter::PrintStyles

Enumerator

***OnlyUUID***

***LOADBULKDATA***

## 27.369.2 Constructor & Destructor Documentation

### 27.369.2.1 gdcm::XMLPrinter::XMLPrinter ( )

### 27.369.2.2 virtual gdcm::XMLPrinter::~~XMLPrinter ( ) [virtual]

## 27.369.3 Member Function Documentation

### 27.369.3.1 PrintStyles gdcm::XMLPrinter::GetPrintStyle ( ) const [inline]

### 27.369.3.2 virtual void gdcm::XMLPrinter::HandleBulkData ( const char \* *uuid*, const TransferSyntax & *ts*, const char \* *bulkdata*, size\_t *bulklen* ) [virtual]

Virtual function mechanism to allow application programmer to override the default mechanism for BulkData handling. By default GDCM will simply discard the BulkData and only write the UUID

### 27.369.3.3 void gdcm::XMLPrinter::Print ( std::ostream & *os* )

### 27.369.3.4 VR gdcm::XMLPrinter::PrintDataElement ( std::ostream & *os*, const Dicts & *dicts*, const DataSet & *ds*, const DataElement & *de*, const TransferSyntax & *ts* ) [protected]

### 27.369.3.5 void gdcm::XMLPrinter::PrintDataSet ( const DataSet & *ds*, const TransferSyntax & *ts*, std::ostream & *os* )

### 27.369.3.6 void gdcm::XMLPrinter::PrintSQ ( const SequenceOfItems \* *sqi*, const TransferSyntax & *ts*, std::ostream & *os* ) [protected]

### 27.369.3.7 void gdcm::XMLPrinter::SetFile ( File const & *f* ) [inline]

### 27.369.3.8 void gdcm::XMLPrinter::SetStyle ( PrintStyles *ps* ) [inline]

## 27.369.4 Member Data Documentation

### 27.369.4.1 const File\* gdcm::XMLPrinter::F [protected]

### 27.369.4.2 PrintStyles gdcm::XMLPrinter::PrintStyle [protected]

The documentation for this class was generated from the following file:


- [gdcmXMLPrinter.h](#)

## 27.370 gdcmm::XMLPrivateDictReader Class Reference


Class for representing a [XMLPrivateDictReader](#).

```
#include <gdcmmXMLPrivateDictReader.h>
```

Inheritance diagram for gdcmm::XMLPrivateDictReader:

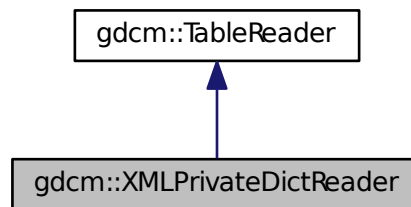

Collaboration diagram for gdcmm::XMLPrivateDictReader:

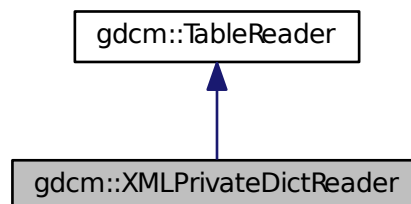

### Public Member Functions

- [XMLPrivateDictReader](#) ()
- [~XMLPrivateDictReader](#) ()
- void [CharacterDataHandler](#) (const char \*data, int length)
- void [EndElement](#) (const char \*name)
- const [PrivateDict](#) & [GetPrivateDict](#) ()
- void [StartElement](#) (const char \*name, const char \*\*atts)

### Protected Member Functions

- void [HandleDescription](#) (const char \*\*atts)
- void [HandleEntry](#) (const char \*\*atts)

### 27.370.1 Detailed Description

Class for representing a [XMLPrivateDictReader](#).

#### Note

bla Will read the Private.xml file

### 27.370.2 Constructor & Destructor Documentation

27.370.2.1 `gdcm::XMLPrivateDictReader::XMLPrivateDictReader ( )`

27.370.2.2 `gdcm::XMLPrivateDictReader::~~XMLPrivateDictReader ( )` `[inline]`

### 27.370.3 Member Function Documentation

27.370.3.1 `void gdcm::XMLPrivateDictReader::CharacterDataHandler ( const char * data, int length )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).

27.370.3.2 `void gdcm::XMLPrivateDictReader::EndElement ( const char * name )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).

27.370.3.3 `const PrivateDict& gdcm::XMLPrivateDictReader::GetPrivateDict ( )` `[inline]`

27.370.3.4 `void gdcm::XMLPrivateDictReader::HandleDescription ( const char ** atts )` `[protected]`

27.370.3.5 `void gdcm::XMLPrivateDictReader::HandleEntry ( const char ** atts )` `[protected]`

27.370.3.6 `void gdcm::XMLPrivateDictReader::StartElement ( const char * name, const char ** atts )` `[virtual]`

Reimplemented from [gdcm::TableReader](#).

The documentation for this class was generated from the following file:

- [gdcmXMLPrivateDictReader.h](#)



## Chapter 28

# File Documentation

### 28.1 gdc2pnm.dox File Reference


### 28.2 gdc2vtk.dox File Reference

### 28.3 gdcmAAbortPDU.h File Reference

```
#include "gdcmTypes.h"
```

```
#include "gdcmBasePDU.h"
```

Include dependency graph for gdcmAAbortPDU.h:

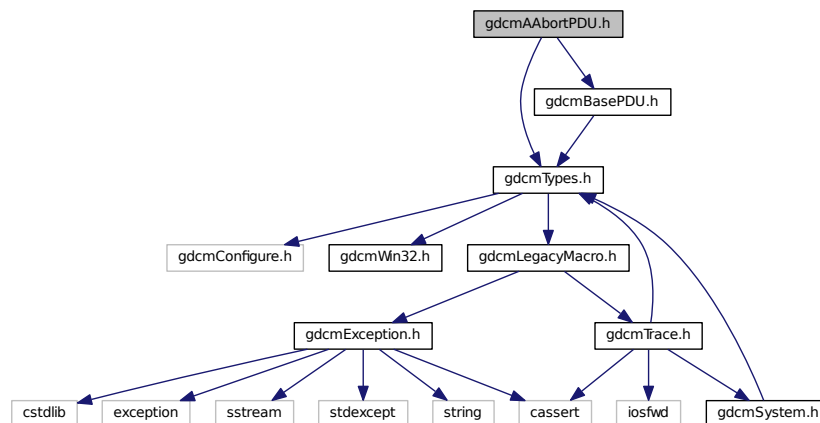

## Classes

- class `gdcm::network::AAbortPDU`

*AAbortPDU Table 9-26 A-ABORT PDU FIELDS.*


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.4 gdcmAAssociateACPDU.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmApplicationContext.h"
#include "gdcmPresentationContextAC.h"
#include "gdcmUserInformation.h"
#include "gdcmBasePDU.h"
#include <vector>
```

Include dependency graph for gdcmAAssociateACPDU.h:

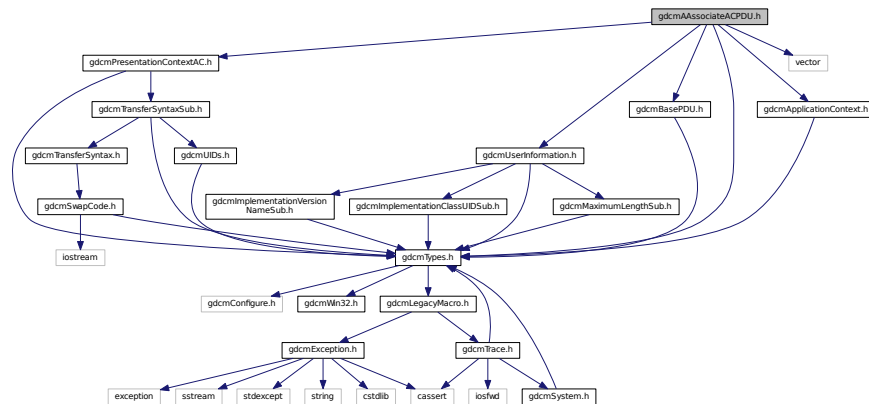

## Classes

- class [gdcm::network::AAssociateACPDU](#)  
*AAssociateACPDU* Table 9-17 ASSOCIATE-AC PDU fields.


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.5 gdcmAAssociateRJPDU.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmBasePDU.h"
```

Include dependency graph for gdcmAAssociateRJPDU.h:

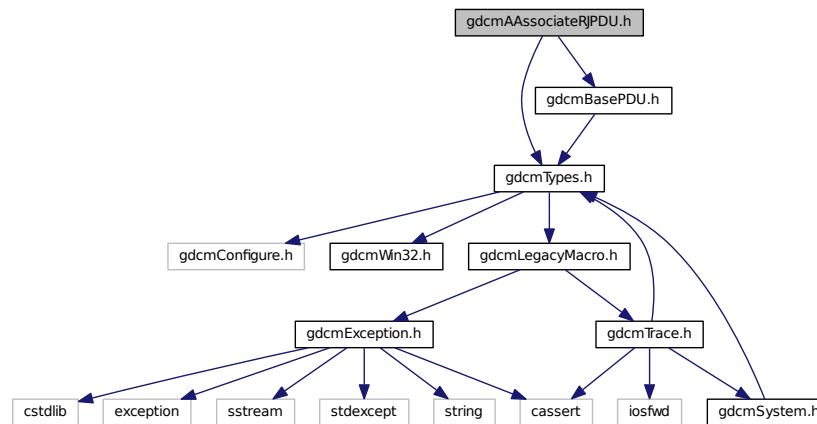

## Classes

- class [gdcm::network::AAssociateRJPDU](#)  
[AAssociateRJPDU](#) Table 9-21 ASSOCIATE-RJ PDU FIELDS.

## Namespaces

- [gdcm](#)
- [gdcm::network](#)


## 28.6 gdcmAAssociateRQPDU.h File Reference

```

#include "gdcmTypes.h"
#include "gdcmVR.h"
#include "gdcmApplicationContext.h"
#include "gdcmPresentationContextRQ.h"
#include "gdcmUserInfo.h"
#include "gdcmBasePDU.h"

```

Include dependency graph for gdcmAAssociateRQPDU.h:

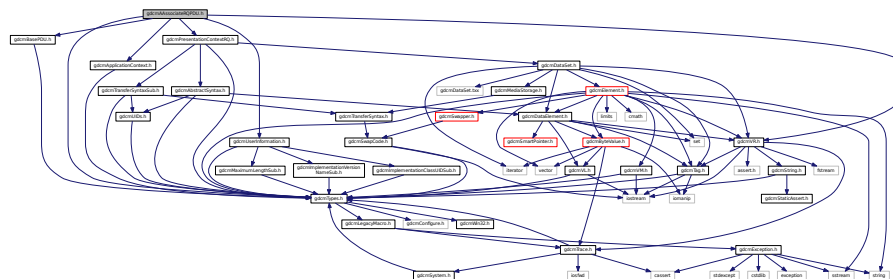

## Classes


- class [gdcmm::network::AAssociateRQPDU](#)  
[AAssociateRQPDU](#) Table 9-11 ASSOCIATE-RQ PDU fields.

## Namespaces


- [gdcmm](#)
- [gdcmm::network](#)

## 28.7 gdcmmAbstractSyntax.h File Reference

```
#include "gdcmmTypes.h"
#include "gdcmmUIDs.h"
#include "gdcmmDataElement.h"
Include dependency graph for gdcmmAbstractSyntax.h:
```

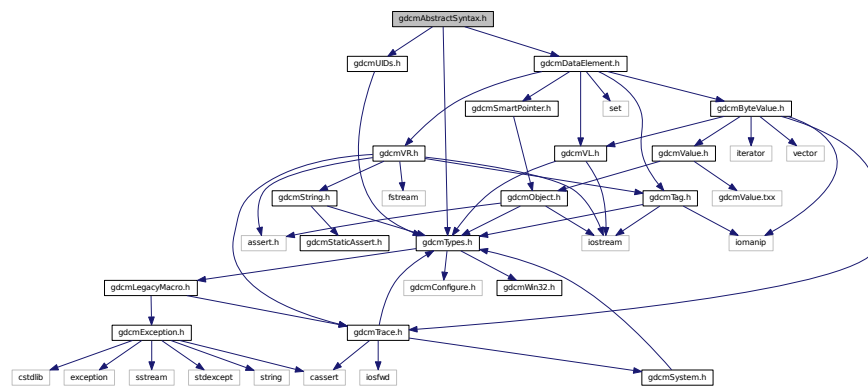

This graph shows which files directly or indirectly include this file:

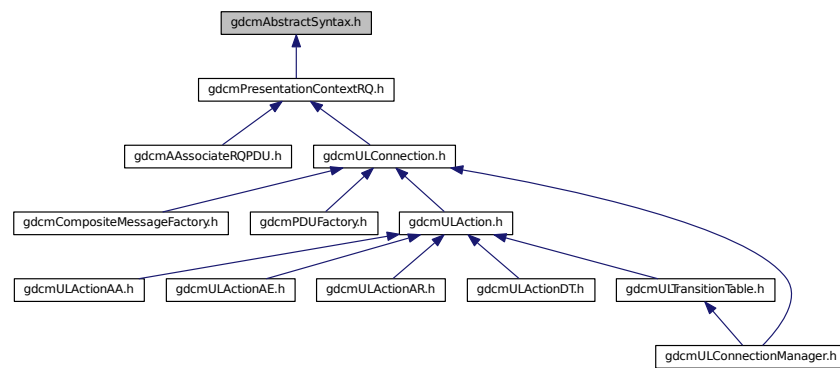



## Classes

- class [gdcm::network::AbstractSyntax](#)

*AbstractSyntax Table 9-14 ABSTRACT SYNTAX SUB-ITEM FIELDS.*

## Namespaces

- [gdcm](#)
- [gdcm::network](#)


## 28.8 gdcmanon.dox File Reference


## 28.9 gdcmAnonymizeEvent.h File Reference

```
#include "gdcmEvent.h"
```

```
#include "gdcmTag.h"
```

Include dependency graph for gdcmAnonymizeEvent.h:

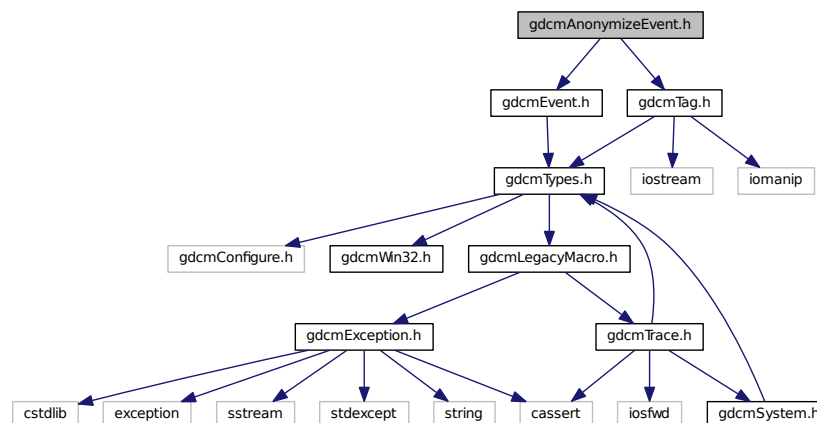

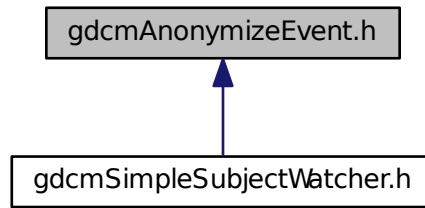


- class `gdcm::AnonymizeEvent`  
*AnonymizeEvent* Special type of event triggered during the Anonymization process.

- **gdcm**

```
#include "gdcmFile.h"
#include "gdcmSubject.h"
#include "gdcmEvent.h"
#include "gdcmSmartPointer.h"
#include <map>
```

Include dependency graph for gdcmAnonymizer.h:

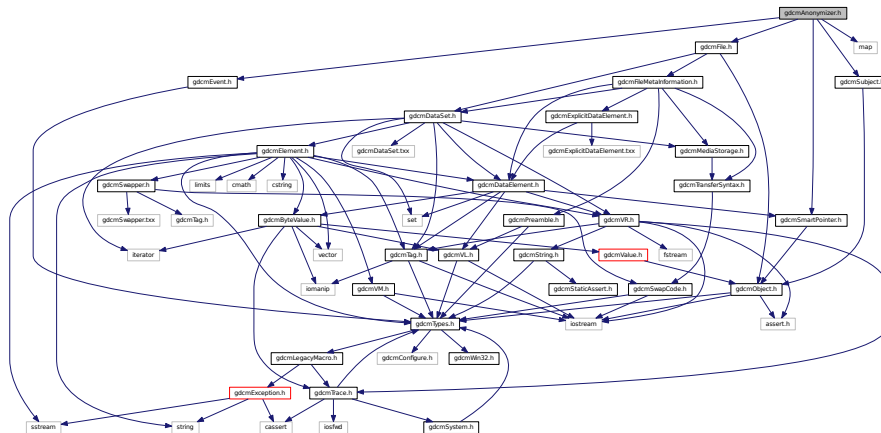

## Classes

- class [gdcm::Anonymizer](#)  
*Anonymizer* This class is a multi purpose anonymizer. It can work in 2 mode:


## Namespaces

- [gdcm](#)


## 28.11 gdcmApplicationContext.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmApplicationContext.h:

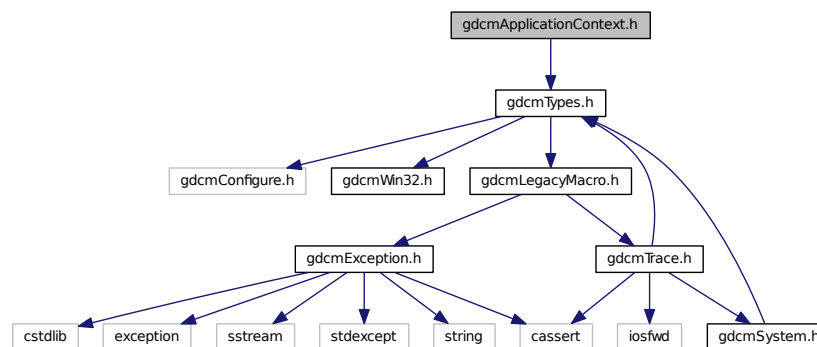

This graph shows which files directly or indirectly include this file:

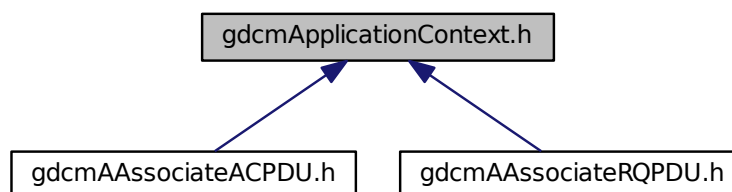

## Classes

- class [gdcm::network::ApplicationContext](#)  
*ApplicationContext* Table 9-12 APPLICATION CONTEXT ITEM FIELDS.


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.12 gdcmApplicationEntity.h File Reference

```
#include "gdcmTypes.h"
#include <vector>
#include <stdlib.h>
```

Include dependency graph for gdcmApplicationEntity.h:

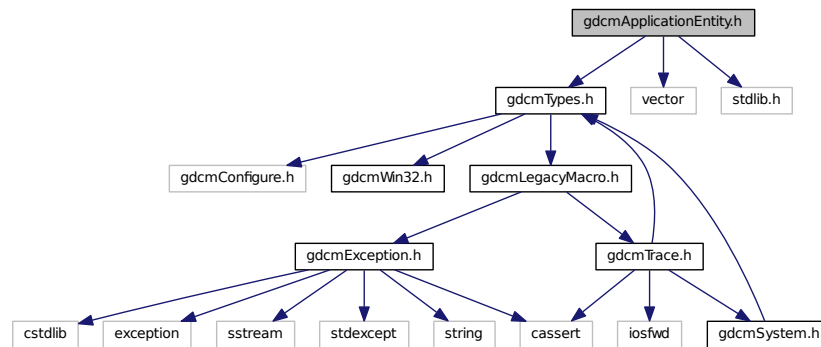

## Classes

- class [gdcm::ApplicationEntity](#)  
*ApplicationEntity.*


## Namespaces

- [gdcm](#)

## 28.13 gdcmAReleaseRPPDU.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmBasePDU.h"
```

Include dependency graph for gdcmAReleaseRPPDU.h:

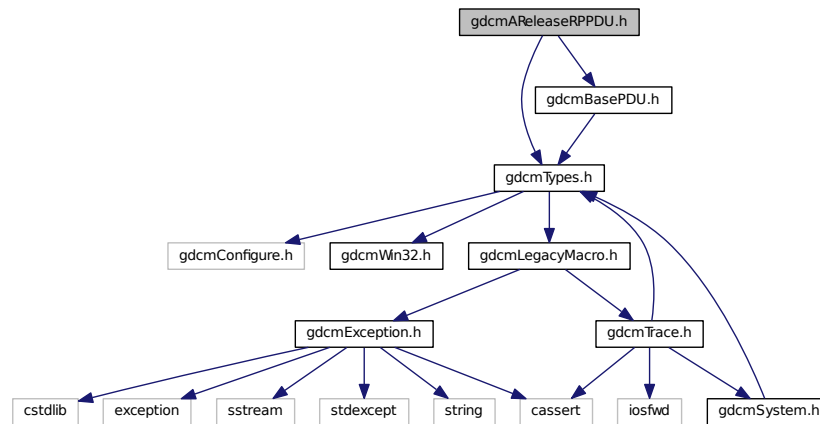

## Classes

- class [gdcm::network::AReleaseRPPDU](#)

[AReleaseRPPDU](#) Table 9-25 A-RELEASE-RP PDU fields.


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.14 gdcmAReleaseRQPDU.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmBasePDU.h"
```

Include dependency graph for `gdcmAReleaseRQPDU.h`:

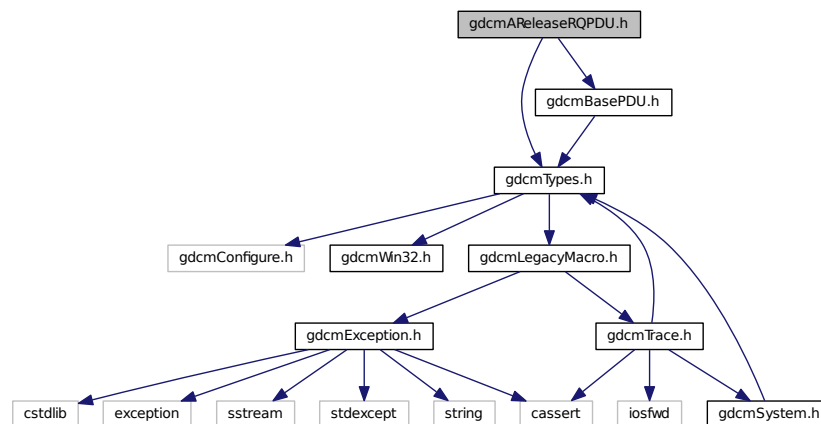

## Classes


- class `gdcm::network::AReleaseRQPDU`  
[AReleaseRQPDU Table 9-24 A-RELEASE-RQ PDU FIELDS.](#)

## Namespaces

- `gdcm`
- `gdcm::network`

## 28.15 gdcmARTIMTimer.h File Reference

This graph shows which files directly or indirectly include this file:

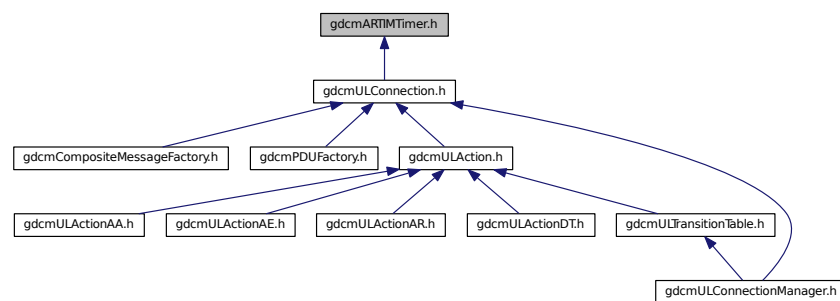

## Classes

- class [gdcm::network::ARTIMTimer](#)  
*ARTIMTimer* This file contains the code for the ARTIM timer.


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.16 gdcmASN1.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmASN1.h:

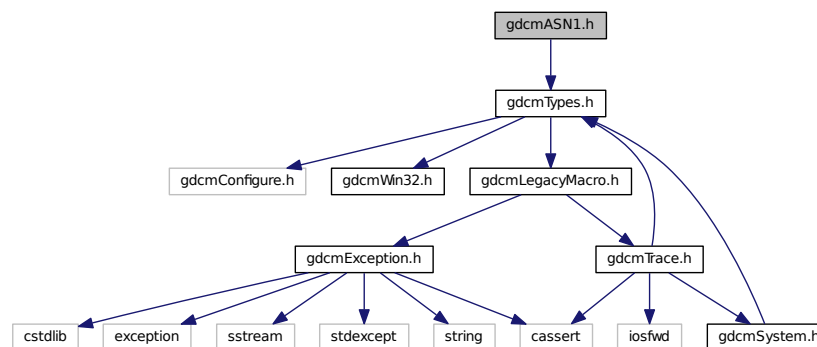

## Classes

- class [gdcm::ASN1](#)  
*Class for ASN1.*


## Namespaces

- [gdcm](#)

## 28.17 gdcmAsynchronousOperationsWindowSub.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmAsynchronousOperationsWindowSub.h`:

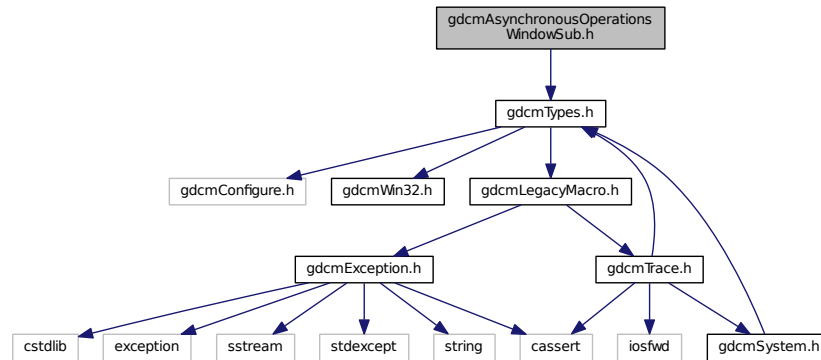

## Classes

- class [gdcm::network::AsynchronousOperationsWindowSub](#)

*AsynchronousOperationsWindowSub* PS 3.7 Table D.3-7 ASYNCHRONOUS OPERATIONS WINDOW SUB-ITEM FIE↔  
LDS (A-ASSOCIATE-RQ)

## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.18 gdcmAttribute.h File Reference


```

#include "gdcmTypes.h"
#include "gdcmVR.h"
#include "gdcmTagToType.h"
#include "gdcmVM.h"
#include "gdcmElement.h"
#include "gdcmDataElement.h"
#include "gdcmDataSet.h"
#include "gdcmStaticAssert.h"
#include <string>
#include <vector>
#include <sstream>


```



Include dependency graph for gdcmAttribute.h:

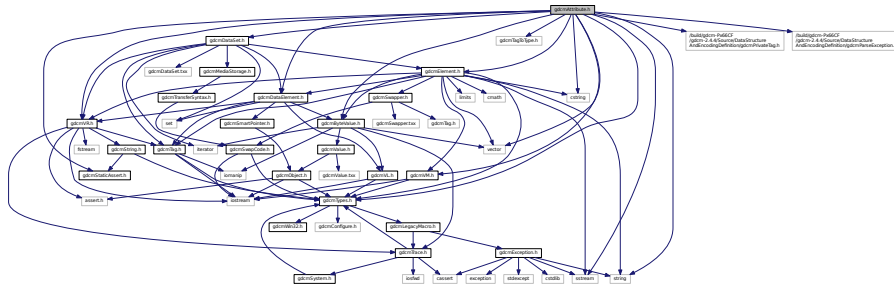

This graph shows which files directly or indirectly include this file:

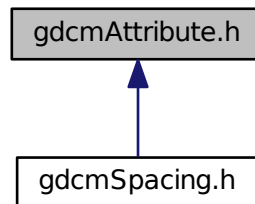

## Classes

- class [gdcm::Attribute< Group, Element, TVR, TVM >](#)  
*Attribute class This class use template metaprograming tricks to let the user know when the template instanciation does not match the public dictionary.*
- class [gdcm::Attribute< Group, Element, TVR, VM::VM1 >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM1\\_3 >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM1\\_8 >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM1\\_n >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM2\\_2n >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM2\\_n >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM3\\_3n >](#)
- class [gdcm::Attribute< Group, Element, TVR, VM::VM3\\_n >](#)
- class [gdcm::VRVLSIZE< T >](#)
- class [gdcm::VRVLSIZE< 0 >](#)
- class [gdcm::VRVLSIZE< 1 >](#)


## Namespaces

- [gdcm](#)

## 28.19 gdcmAudioCodec.h File Reference

```
#include "gdcmCodec.h"
```

Include dependency graph for gdcmAudioCodec.h:

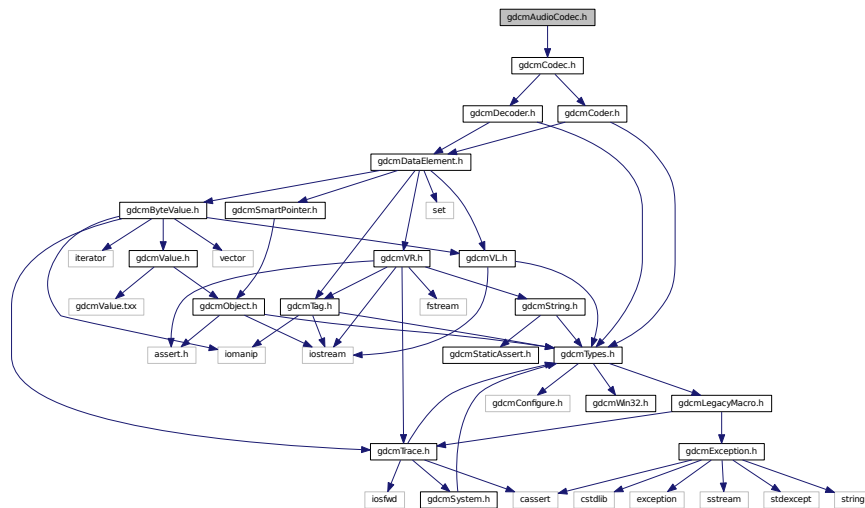

### Classes

- class [gdcm::AudioCodec](#)

*AudioCodec.*


### Namespaces

- [gdcm](#)

## 28.20 gdcmBase64.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmbase64.h:

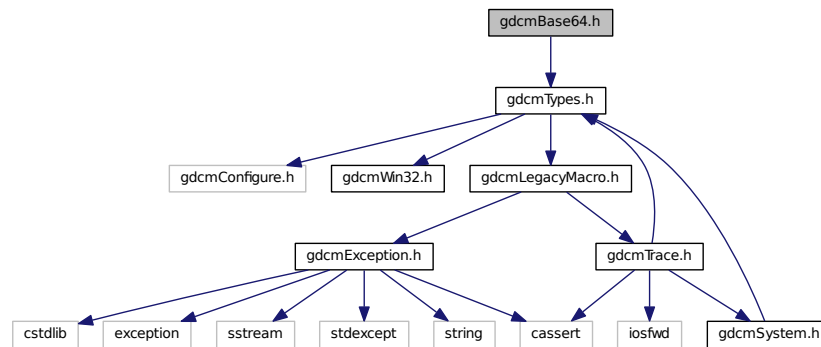

## Classes

- class `gdcm::Base64`  
*Class for Base64.*


## Namespaces

- **gdcm**


## 28.21 gdcmBaseCompositeMessage.h File Reference

```
#include "gdcmPresentationDataValue.h"
#include "gdcmBaseRootQuery.h"
#include <vector>
```

Include dependency graph for `gdcmBaseCompositeMessage.h`:

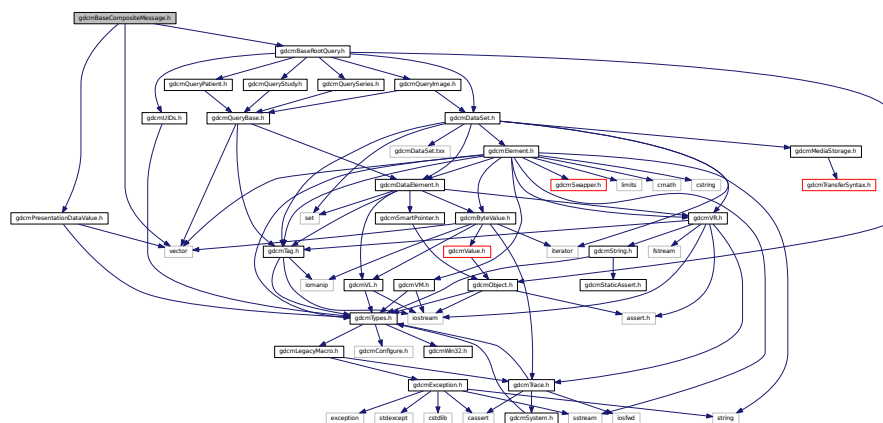

This graph shows which files directly or indirectly include this file:

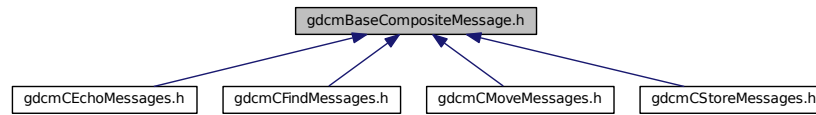

## Classes

- class [gdcmbaseCompositeMessage](#)

*[BaseCompositeMessage](#) The Composite events described in section 3.7-2009 of the DICOM standard all use their own messages. These messages are constructed using Presentation Data Values, from section 3.8-2009 of the standard, and then fill in appropriate values in their datasets.*


## Namespaces

- [gdcmbase](#)
- [gdcmbase::network](#)

## 28.22 gdcmbasePDU.h File Reference

```
#include "gdcmbaseTypes.h"
```

Include dependency graph for `gdcmbasePDU.h`:

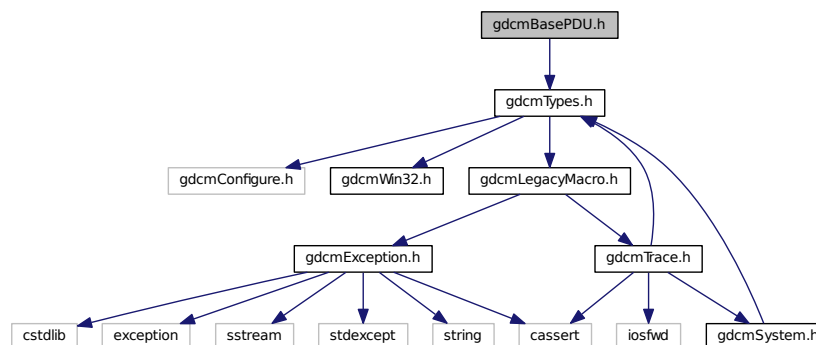


```

graph TD
    gbmmlibPOU[gbmmlib.POU.h] --> gbmmlibActorPOU[gbmmlibActor.POU.h]
    gbmmlibPOU --> gbmmlibAccumulatorPOU[gbmmlibAccumulator.POU.h]
    gbmmlibPOU --> gbmmlibAccumulator2POU[gbmmlibAccumulator2.POU.h]
    gbmmlibPOU --> gbmmlibAccumulator3POU[gbmmlibAccumulator3.POU.h]
    gbmmlibPOU --> gbmmlibAccumulator4POU[gbmmlibAccumulator4.POU.h]
    gbmmlibPOU --> gbmmlibAccumulator5POU[gbmmlibAccumulator5.POU.h]
    gbmmlibPOU --> gbmmlibDataPOU[gbmmlibData.POU.h]
    gbmmlibPOU --> gbmmlibEvent[gbmmlibEvent.h]
    gbmmlibEvent --> gbmmlibAction[gbmmlibAction.h]
    gbmmlibAction --> gbmmlibActionAA[gbmmlibActionAA.h]
    gbmmlibAction --> gbmmlibActionAE[gbmmlibActionAE.h]
    gbmmlibAction --> gbmmlibActionAB[gbmmlibActionAB.h]
    gbmmlibAction --> gbmmlibActionDTH[gbmmlibActionDTH.h]
    gbmmlibAction --> gbmmlibTransitionTable[gbmmlibTransitionTable.h]
    gbmmlibTransitionTable --> gbmmlibConnectorManager[gbmmlibConnectorManager.h]
  
```


- class `gdcn::network::BasePDU`  
*BasePDU* base class for PDUs.

- `gdcm`
- `gdcm::network`

```
#include "gdcmDataSet.h"
#include "gdcmUIDs.h"
#include "gdcmObject.h"
#include "gdcmQueryPatient.h"
#include "gdcmQueryStudy.h"
#include "gdcmQuerySeries.h"
#include "gdcmQueryImage.h"
Include dependency graph for gdcmBaseRootQuery.h:
```

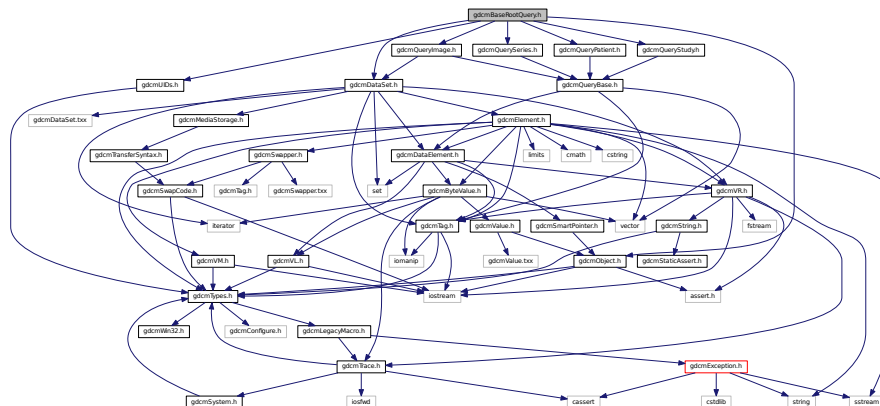

This graph shows which files directly or indirectly include this file:

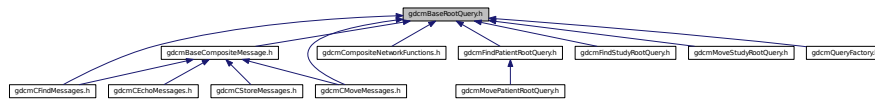

## Classes

- class [gdcmbase::BaseRootQuery](#)

*BaseRootQuery* contains: a baseclass which will produce a dataset for c-find and c-move with patient/study root.

## Namespaces

- [gdcmbase](#)


## Enumerations

- enum [gdcmbase::EQueryLevel](#) {  
[gdcmbase::ePatient](#) = 0,  
[gdcmbase::eStudy](#) = 1,  
[gdcmbase::eSeries](#) = 2,  
[gdcmbase::eImage](#) = 3 }
- enum [gdcmbase::EQueryType](#) {  
[gdcmbase::eFind](#) = 0,  
[gdcmbase::eMove](#) }


## 28.24 gdcmbasicOffsetTable.h File Reference

```
#include "gdcmbaseFragment.h"
```

Include dependency graph for gdcmBasicOffsetTable.h:

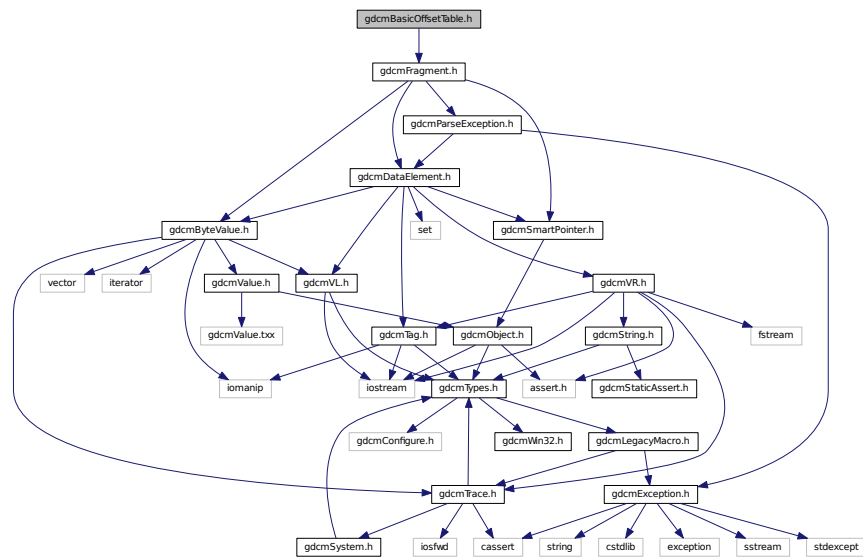

This graph shows which files directly or indirectly include this file:

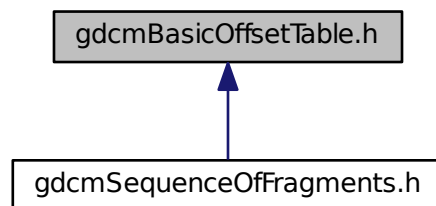

## Classes

- class [gdcm::BasicOffsetTable](#)  
Class to represent a *BasicOffsetTable*.

## Namespaces

- [gdcm](#)






## Namespaces

- **gdcm**


## 28.26 gdcmBitmapToBitmapFilter.h File Reference

```
#include "gdcmBitmap.h"
```

Include dependency graph for gdcmBitmapToBitmapFilter.h:

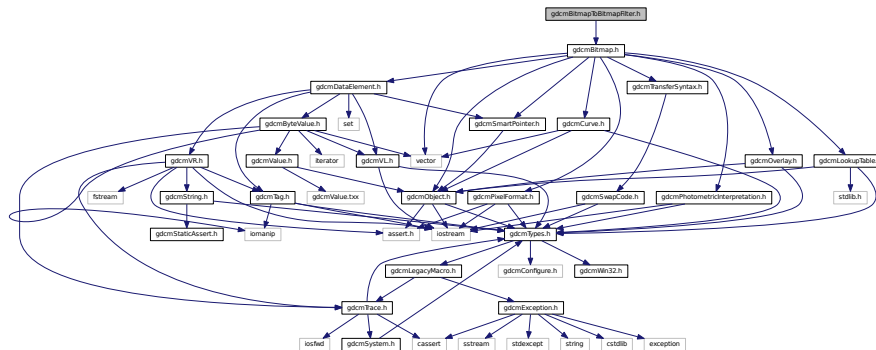

This graph shows which files directly or indirectly include this file:

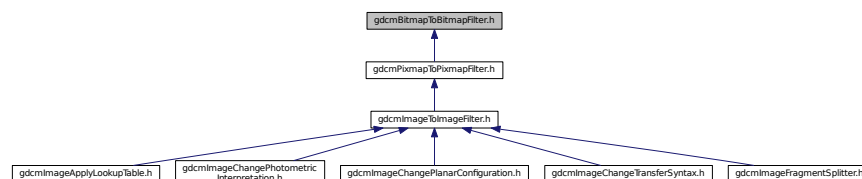

## Classes

- class `gdcm::BitmapToBitmapFilter`  
*BitmapToBitmapFilter* class Super class for all filter taking an image and producing an output image.


## Namespaces

- gdc

## 28.27 gdcmBoxRegion.h File Reference

```
#include "gdcmRegion.h"
```

Include dependency graph for `gdcmBoxRegion.h`:

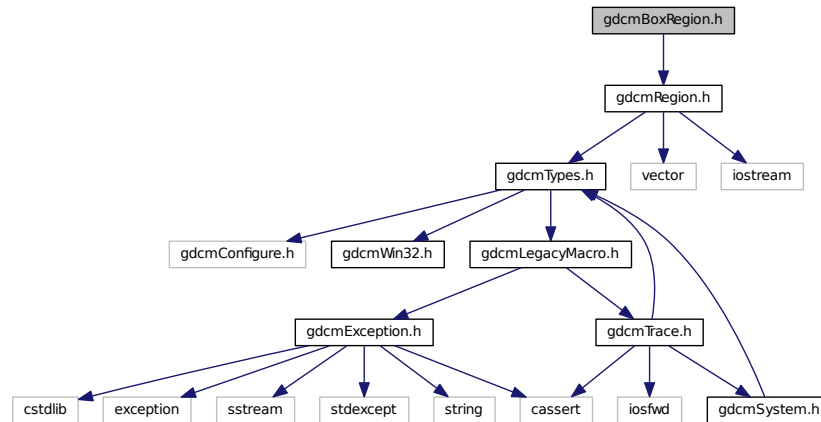

## Classes

- class [gdcm::BoxRegion](#)

*Class for manipulation box region This is a very simple implementation of the [Region](#) class. It only support 3D box type region. It assumes the 3D Box does not have a tilt Origin is as (0,0,0)*

## Namespaces

- [gdcm](#)


## 28.28 gdcmByteBuffer.h File Reference

```


#include "gdcmTypes.h"
#include <vector>
#include <assert.h>
#include <string.h>

```

Include dependency graph for gdcmByteBuffer.h:

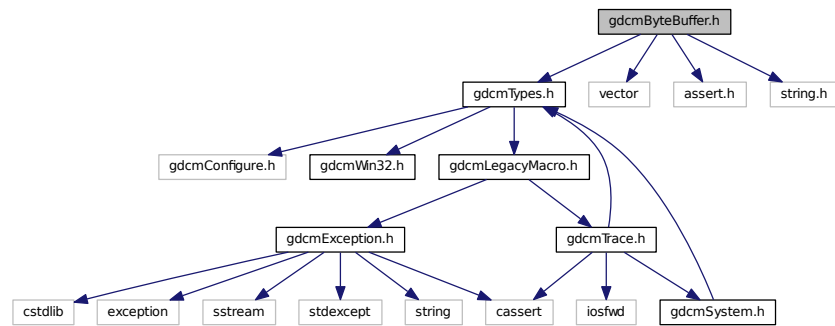

This graph shows which files directly or indirectly include this file:

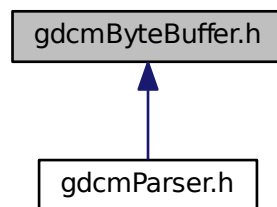

## Classes

- class [gdcm::ByteBuffer](#)  
*ByteBuffer.*

## Namespaces


- [gdcm](#)

## 28.29 gdcmByteSwap.h File Reference

```

#include "gdcmTypes.h"
#include "gdcmSwapCode.h"
#include "gdcmByteSwap.txx"


```

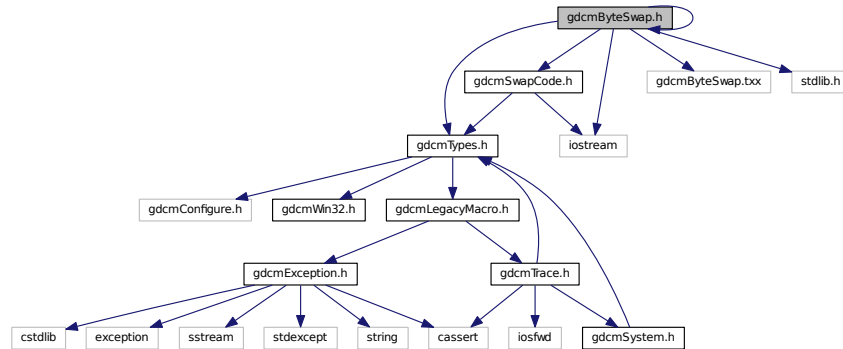

- class `gdcm::ByteSwap< T >`  
*ByteSwap.*

- **gdcm**

```
#include "gdcmDataSet.h"
```

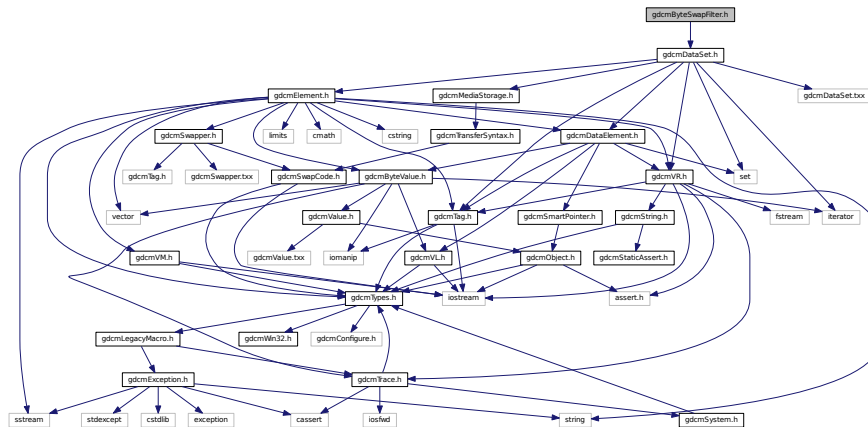

## Classes

- class [gdcm::ByteSwapFilter](#)  
*ByteSwapFilter* In place byte-swapping of a dataset *FIXME: FL status ??*


## Namespaces

- [gdcm](#)


## 28.31 gdcmByteValue.h File Reference

```
#include "gdcmValue.h"
#include "gdcmTrace.h"
#include "gdcmVL.h"
#include <vector>
#include <iterator>
#include <iomanip>
```

Include dependency graph for gdcmByteValue.h:

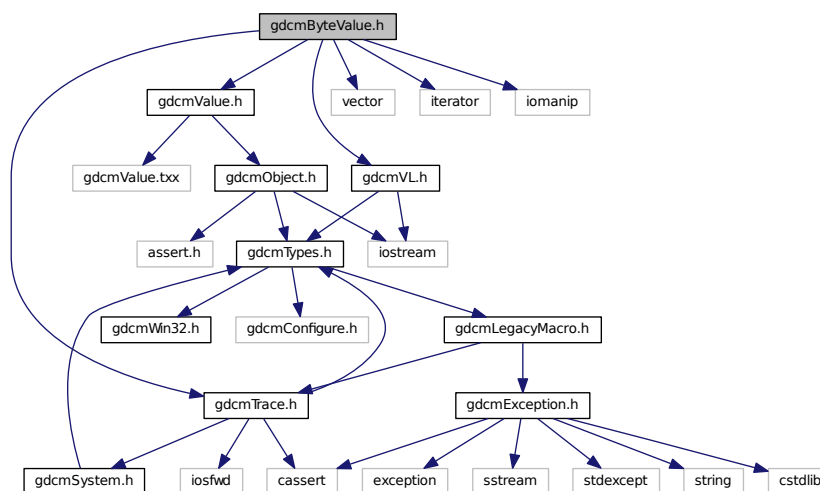

This graph shows which files directly or indirectly include this file:

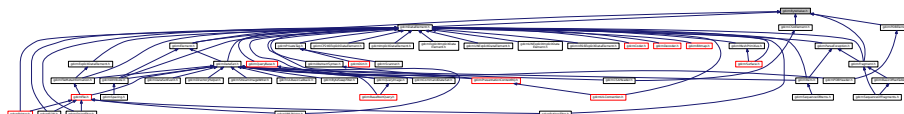

## Classes

- class [gdcm::ByteValue](#)

*Class to represent binary value (array of bytes)*


## Namespaces

- [gdc](#)

## 28.32 gdcmCAPICryptoFactory.h File Reference

```
#include "gdcmCryptoFactory.h"
```

Include dependency graph for gdcmCAPICryptoFactory.h:

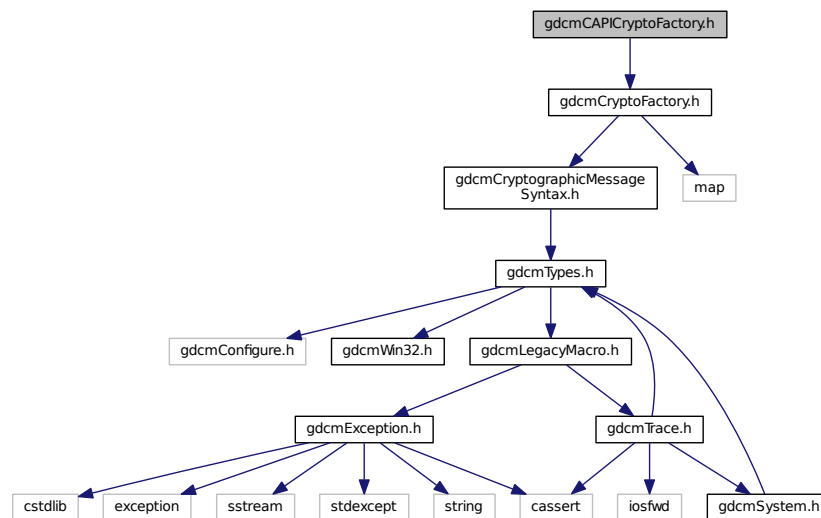

## Classes

- class [gdc::CAPICryptoFactory](#)

## Namespaces

- [gdc](#)

## 28.33 gdcmCAPICryptographicMessageSyntax.h File Reference

```
#include "gdcmCryptographicMessageSyntax.h"
#include <windows.h>
#include <wincrypt.h>
#include <vector>
```


```

graph TD
    Root["gdcmCAPICryptographicMessageSyntax.h"] --> GDCMMsgSyntax["gdcmCryptographicMessageSyntax.h"]
    Root --> Windows["windows.h"]
    Root --> Wincrypt["wincrypt.h"]
    Root --> Vector["vector"]
    GDCMMsgSyntax --> GDCMTypes["gdcmTypes.h"]
    GDCMTypes --> GDCMConfigure["gdcmConfigure.h"]
    GDCMTypes --> GDCMWin32["gdcmWin32.h"]
    GDCMTypes --> GDCMLegacyMacro["gdcmLegacyMacro.h"]
    GDCMLegacyMacro --> GDCMException["gdcmException.h"]
    GDCMLegacyMacro --> GDCMTrace["gdcmTrace.h"]
    GDCMException --> Cstdlib["cstdlib"]
    GDCMException --> Exception["exception"]
    GDCMException --> Sstream["sstream"]
    GDCMException --> Stdexcept["stdexcept"]
    GDCMException --> String["string"]
    GDCMException --> Cassert["cassert"]
    GDCMTrace --> IOSfwd["iosfwd"]
    GDCMTrace --> GDCMSystem["gdcmSystem.h"]
    GDCMTypes --> GDCMSystem
  
```

- class `gdc::CAPICryptographicMessageSyntax`

- `gdcm`

```
#include "gdcmBaseCompositeMessage.h"
Include dependency graph for gdcmCEchoMessages.h:
```

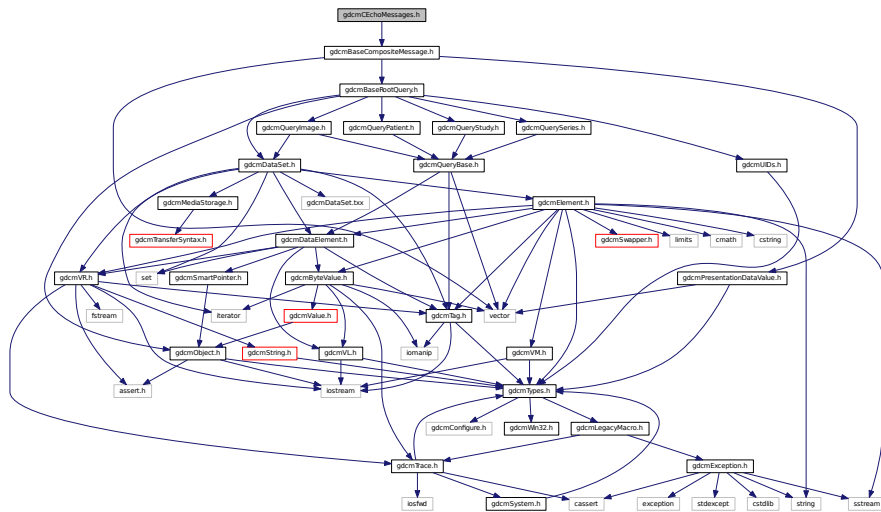









Include dependency graph for `gdcmCodec.h`:

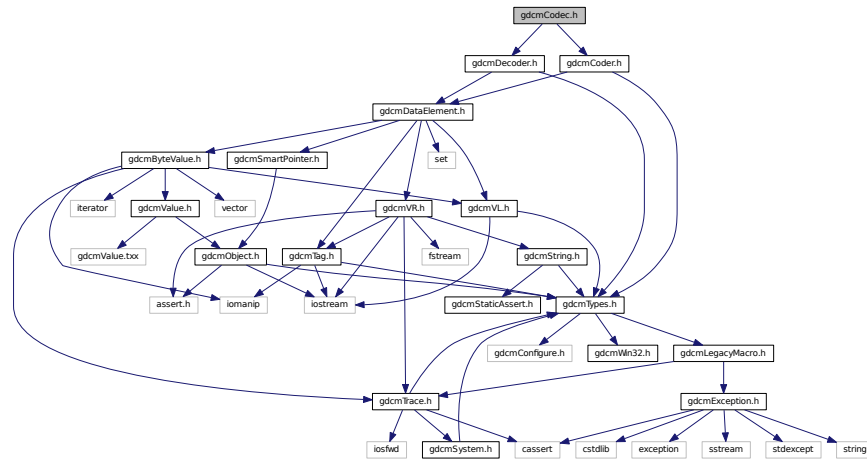

This graph shows which files directly or indirectly include this file:

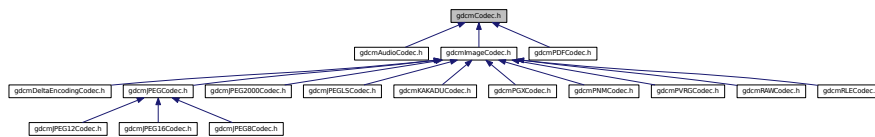

## Classes

- class `gdcm::Codec`

*Codec* class.


## Namespaces

- `gdcm`


## 28.38 gdcmCoder.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmDataElement.h"
```

Include dependency graph for gdcmCoder.h:

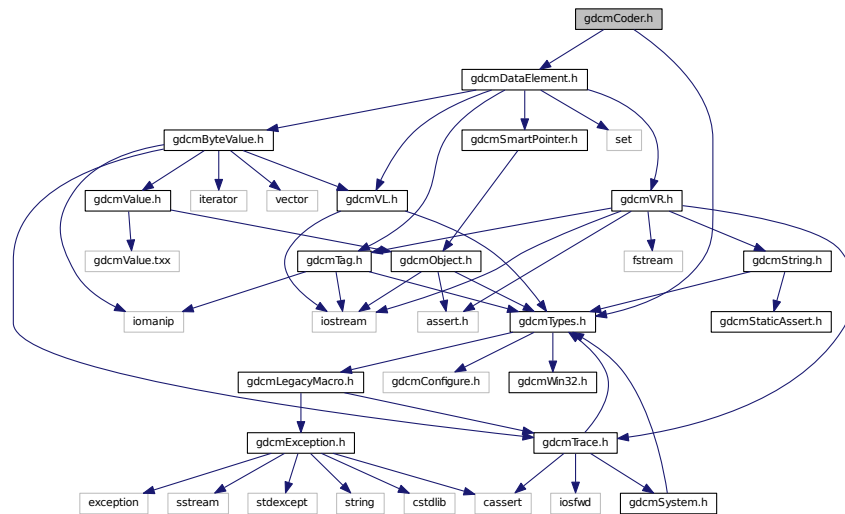

This graph shows which files directly or indirectly include this file:

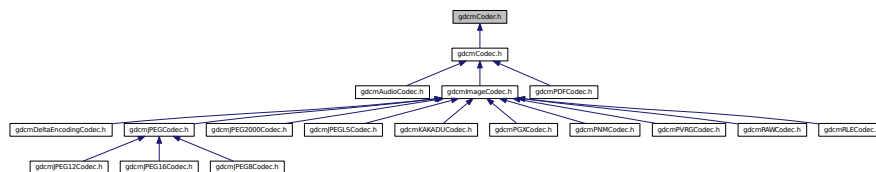

## Classes

- class [gdcm::Coder](#)  
*Coder.*


## Namespaces

- [gdcm](#)

## 28.39 gdcmCodeString.h File Reference

```
#include "gdcmString.h"
```

Include dependency graph for `gdcmCodeString.h`:

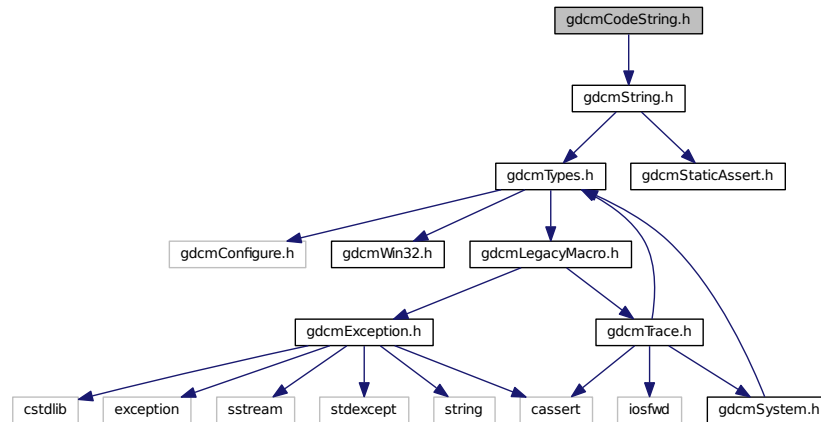

## Classes

- class [gdcm::CodeString](#)

*[CodeString](#) This is an implementation of DICOM [VR](#): CS The ctor will properly Trim so that operator== is correct.*

## Namespaces

- [gdcm](#)


## Functions

- bool [gdcm::operator!=](#) (const CodeString &ref, const CodeString &cs)
- std::ostream & [gdcm::operator<<](#) (std::ostream &os, const CodeString &str)
- bool [gdcm::operator==](#) (const CodeString &ref, const CodeString &cs)


## 28.40 gdcmCommand.h File Reference

```
#include "gdcmSubject.h"
```

Include dependency graph for gdcMCommand.h:

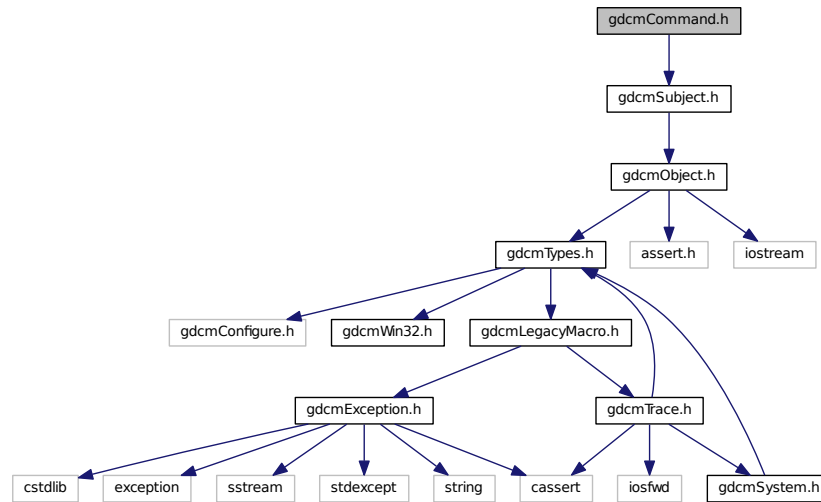

This graph shows which files directly or indirectly include this file:

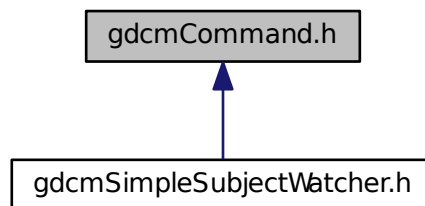

## Classes

- class [gdcM::Command](#)  
*Command superclass for callback/observer methods.*
- class [gdcM::MemberCommand< T >](#)  
*Command subclass that calls a pointer to a member function.*
- class [gdcM::SimpleMemberCommand< T >](#)  
*Command subclass that calls a pointer to a member function.*

## Namespaces


- [gdcM](#)

## 28.41 gdcmCommandDataSet.h File Reference

```
#include "gdcmDataSet.h"
```

```
#include "gdcmDataElement.h"
```

Include dependency graph for gdcmCommandDataSet.h:

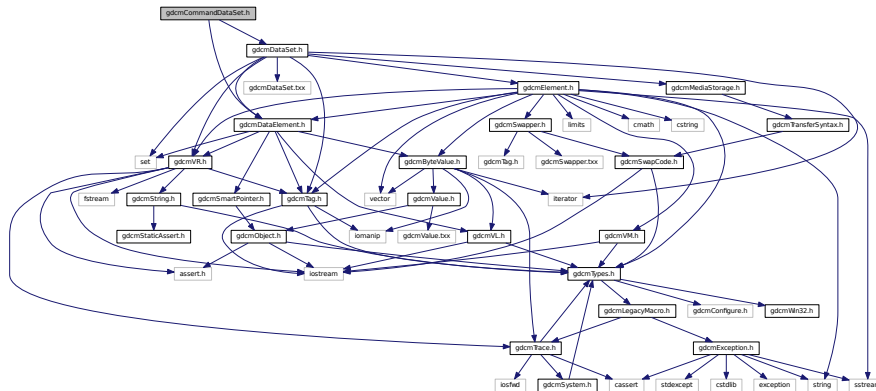

### Classes

- class [gdcm::CommandDataSet](#)

*Class to represent a [Command DataSet](#).*

### Namespaces

- [gdcm](#)

### Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const CommandDataSet &val)`

## 28.42 gdcmCompositeMessageFactory.h File Reference

```
#include "gdcmPresentationDataValue.h"
```

```
#include "gdcmULConnection.h"
```

- class `gdcm::network::CompositeMessageFactory`

## Namespaces

- ## 28.43 gdcmCompositeNetworkFunctions.h File Reference

Generated on Sat Sep 26 2015 07:03:31 for GDCM by Doxygen





## 28.46 gdcmCP246ExplicitDataElement.h File Reference

```
#include "gdcmDataElement.h"
#include "gdcmCP246ExplicitDataElement.txx"
Include dependency graph for gdcmCP246ExplicitDataElement.h:
```



### Classes


- class [gdcm::CP246ExplicitDataElement](#)  
Class to read/write a [DataElement](#) as CP246Explicit Data [Element](#).

### Namespaces


- [gdcm](#)

## 28.47 gdcmCryptoFactory.h File Reference

```
#include "gdcmCryptographicMessageSyntax.h"
#include <map>
Include dependency graph for gdcmCryptoFactory.h:
```

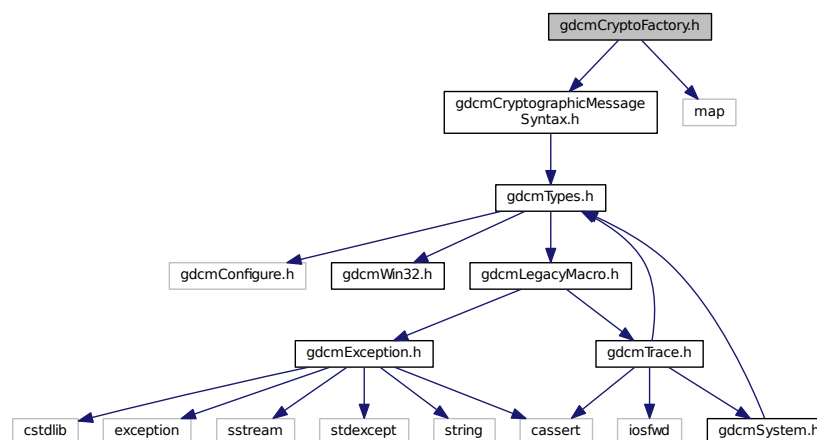

This graph shows which files directly or indirectly include this file:

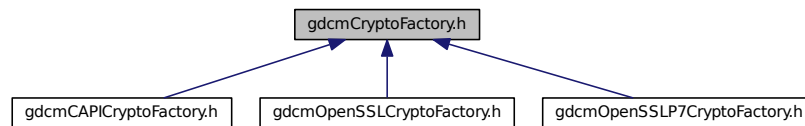

## Classes

- class `gdcm::CryptoFactory`

*Class to do handle the crypto factory.*


## Namespaces

- `gdcm`


## 28.48 gdcmCryptographicMessageSyntax.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmCryptographicMessageSyntax.h`:

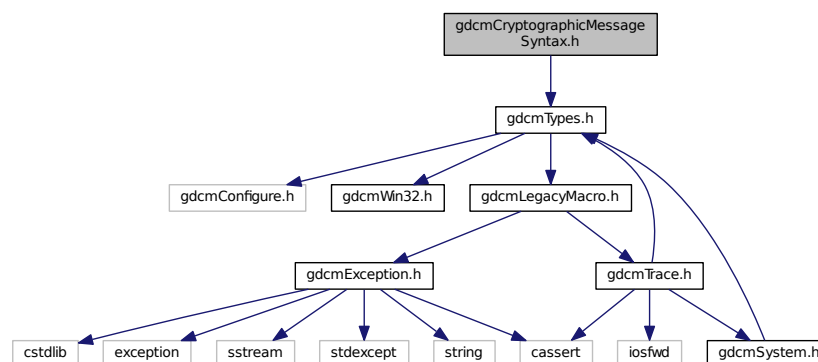

This graph shows which files directly or indirectly include this file:

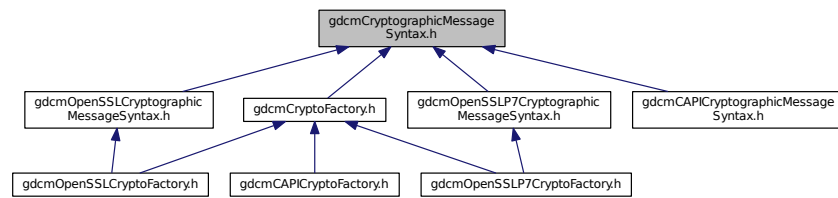

## Classes


- class [gdcm::CryptographicMessageSyntax](#)

## Namespaces


- [gdcm](#)

## 28.49 gdcmCSAElement.h File Reference

```
#include "gdcmTag.h"
#include "gdcmVM.h"
#include "gdcmVR.h"
#include "gdcmByteValue.h"
#include "gdcmSmartPointer.h"
Include dependency graph for gdcmCSAElement.h:
```

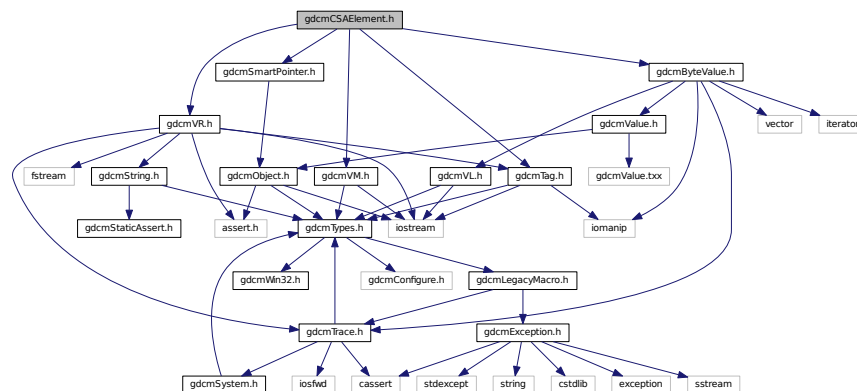

This graph shows which files directly or indirectly include this file:

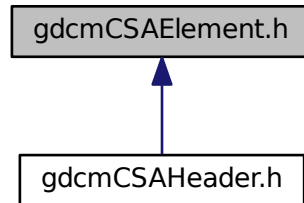

## Classes

- class [gdcm::CSAElement](#)

*Class to represent a CSA [Element](#).*

## Namespaces

- [gdcm](#)

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const CSAElement &val)`

## 28.50 gdcmCSAHeader.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmDataSet.h"
#include "gdcmCSAElement.h"
```

[illegible]

- class `gdcm::CSAHeader`


## Namespaces

- **gdcm**


- `std::ostream & gdcmm::operator<< (std::ostream &os, const CSAHeader &d)`

```
#include "gdcmTypes.h"
#include "gdcmTag.h"
#include "gdcmCSAHeaderDictEntry.h"
#include <iostream>
#include <iomanip>
#include <set>
#include <exception>
```

Include dependency graph for `gdcmsAHeaderDict.h`:

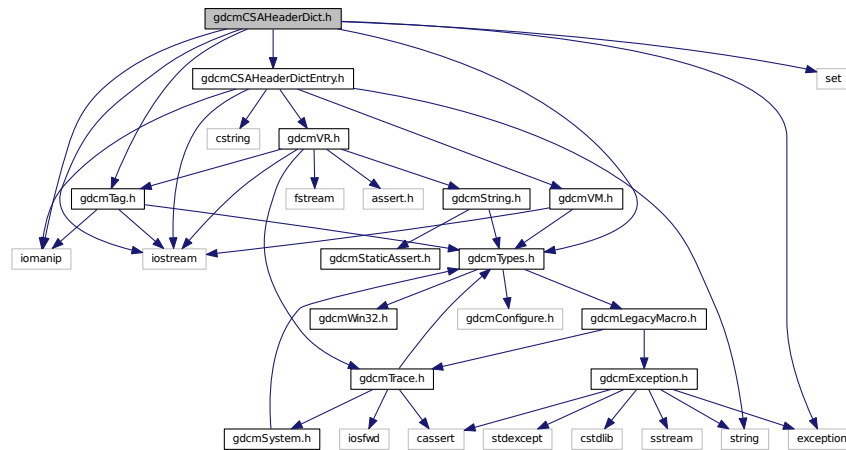

This graph shows which files directly or indirectly include this file:

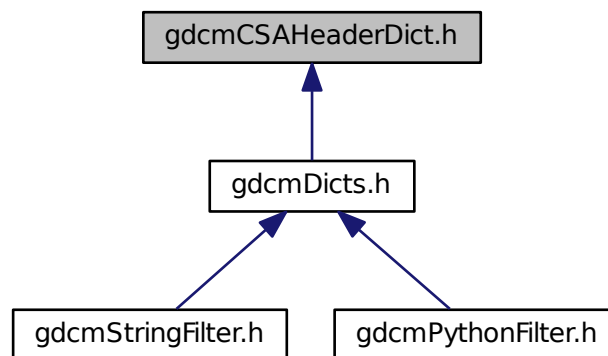

## Classes

- class [gdcms::CSAHeaderDict](#)  
Class to represent a map of [CSAHeaderDictEntry](#).
- class [gdcms::CSAHeaderDictException](#)

## Namespaces

- [gdcms](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const CSAHeaderDict &val)`


## 28.52 gdcmCSAHeaderDictEntry.h File Reference

```
#include "gdcmVR.h"
#include "gdcmVM.h"
#include <string>
#include <iostream>
#include <iomanip>
#include <cstring>
```

Include dependency graph for gdcmCSAHeaderDictEntry.h:

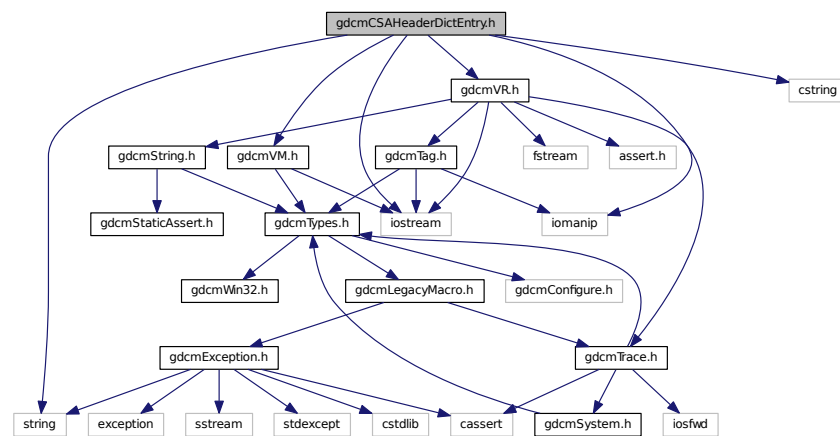

This graph shows which files directly or indirectly include this file:

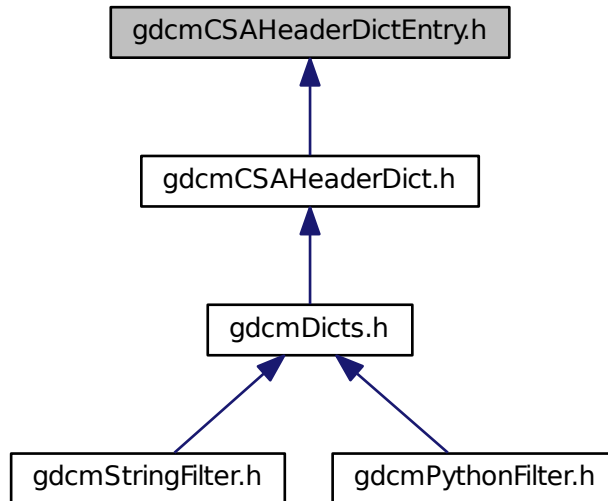

## Classes

- class [gdcml::CSAHeaderDictEntry](#)

*Class to represent an Entry in the [Dict](#) Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcml::Tag](#) to the needed information.*

## Namespaces

- [gdcml](#)

## Functions

- `std::ostream & gdcml::operator<< (std::ostream &os, const CSAHeaderDictEntry &val)`

## 28.53 gdcmlCStoreMessages.h File Reference



```
#include "gdcmlBaseCompositeMessage.h"
```



- class `gdcn::network::CStoreRQ`  
*CStoreRQ* this file defines the messages for the cecho action.
- class `gdcn::network::CStoreRSP`  
*CStoreRSP* this file defines the messages for the cecho action.

- `gdcm`
- `gdcm::network`


```
#include "gdcmTypes.h"
#include "gdcmObject.h"
#include <vector>
```

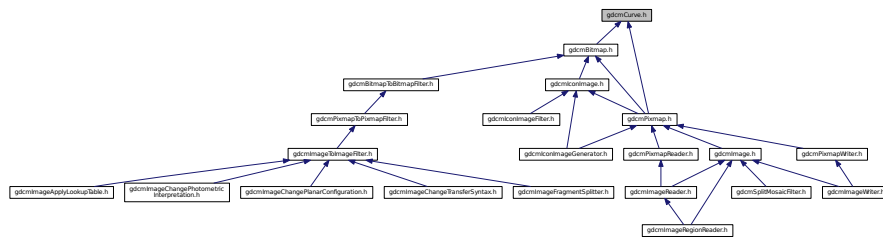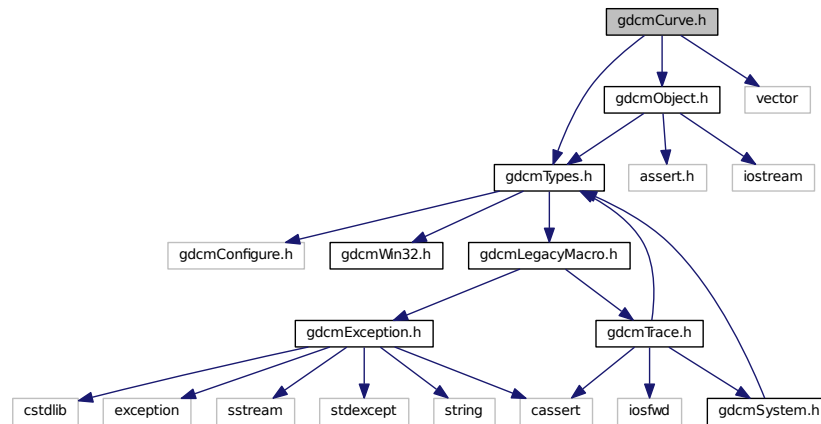


- Curve class to ha

```
#include "gdcmTag.h"
#include "gdcmVL.h"
#include "gdcmVR.h"
#include "gdcmByteValue.h"
#include "gdcmSmartPointer.h"
#include <set>
```

Include dependency graph for gdcmDataElement.h:

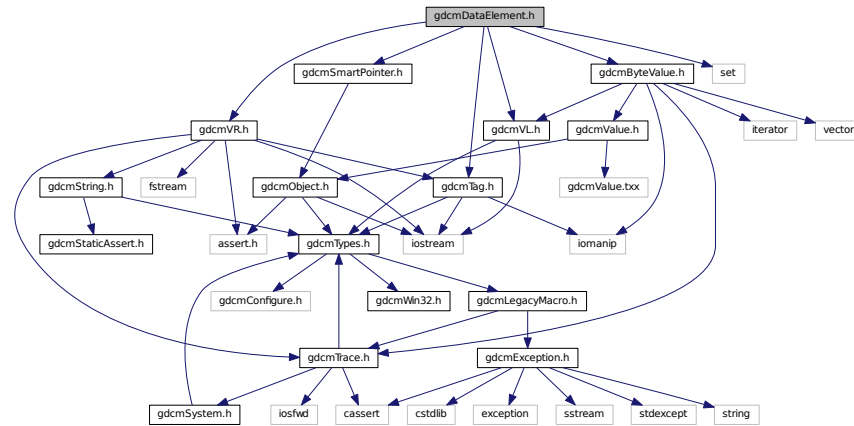

This graph shows which files directly or indirectly include this file:

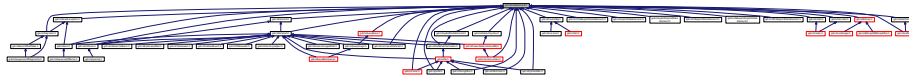

## Classes

- class [gdcm::DataElement](#)  
Class to represent a Data [Element](#) either Implicit or Explicit.

## Namespaces

- [gdcm](#)


## Functions

- bool [gdcm::operator!=](#) (const DataElement &lhs, const DataElement &rhs)
- std::ostream & [gdcm::operator<<](#) (std::ostream &os, const DataElement &val)


## 28.56 gdcmDataEvent.h File Reference

```
#include "gdcmEvent.h"
```

Include dependency graph for `gdcmDataEvent.h`:

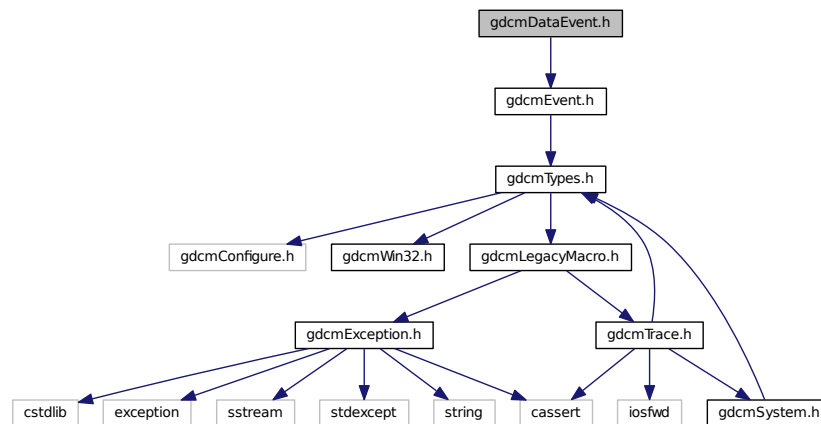

This graph shows which files directly or indirectly include this file:

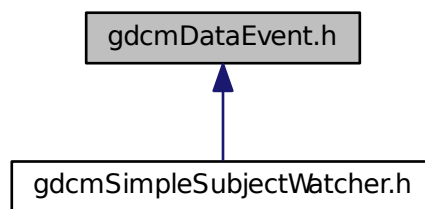

## Classes

- class `gdcm::DataEvent`  
*DataEvent.*


## Namespaces

- `gdcm`


## 28.57 gdcmDataSet.h File Reference

```
#include "gdcmDataElement.h"
```

```
#include "gdcmTag.h"
#include "gdcmVR.h"
#include "gdcmElement.h"
#include "gdcmMediaStorage.h"
#include <set>
#include <iterator>
#include "gdcmDataSet.txx"
Include dependency graph for gdcmDataSet.h:
```

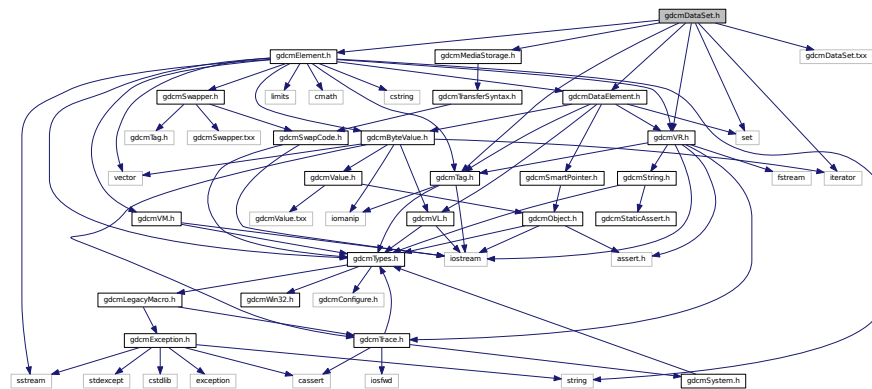

This graph shows which files directly or indirectly include this file:

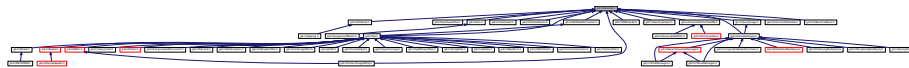

## Classes

- class `gdcm::DataElementException`
- class `gdcm::DataSet`

*Class to represent a Data Set (which contains Data Elements) A Data Set represents an instance of a real world Information Object.*

## Namespaces

- **gdcm**


## Functions

- `std::ostream & gdcmm::operator<< (std::ostream &os, const DataSet &val)`

## 28.58 gdcmDataSetEvent.h File Reference

```
#include "gdcmEvent.h"
```

Include dependency graph for `gdcmDataSetEvent.h`:

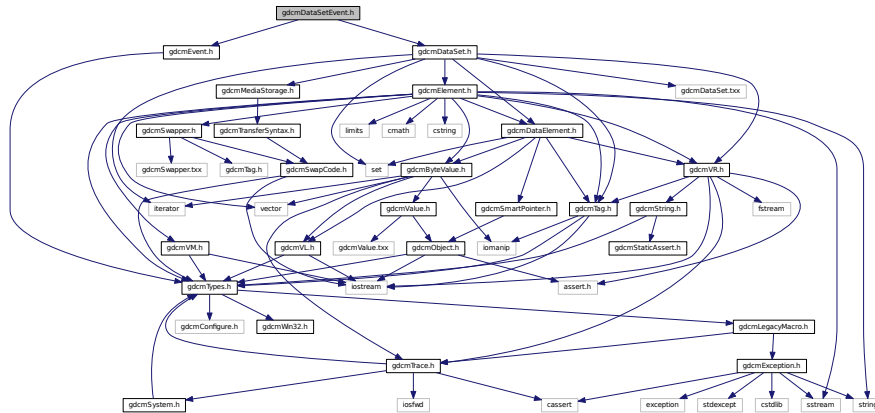

## Classes

- class `gdcm::DataSetEvent`

**DataSetEvent** Special type of event triggered during the **DataSet** store/move process.


## Namespaces

- `gdcm`

## 28.59 gdcmDataSetHelper.h File Reference

```
#include "gdcMVR.h"
```

Include dependency graph for gdcmDataSetHelper.h:

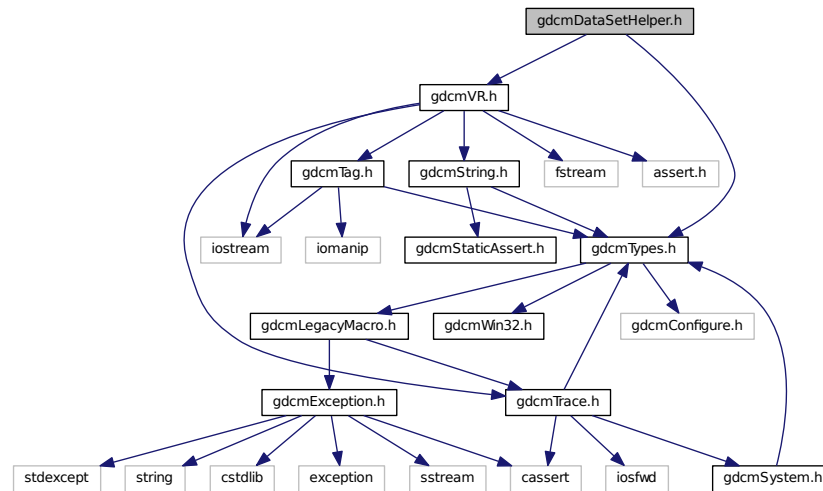

## Classes

- class `gdcm::DataSetHelper`

*DataSetHelper* (internal class, not intended for user level)


## Namespaces

- `gdcm`


## 28.60 gdcmDecoder.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmDataElement.h"
```

Include dependency graph for `gdcmDecoder.h`:

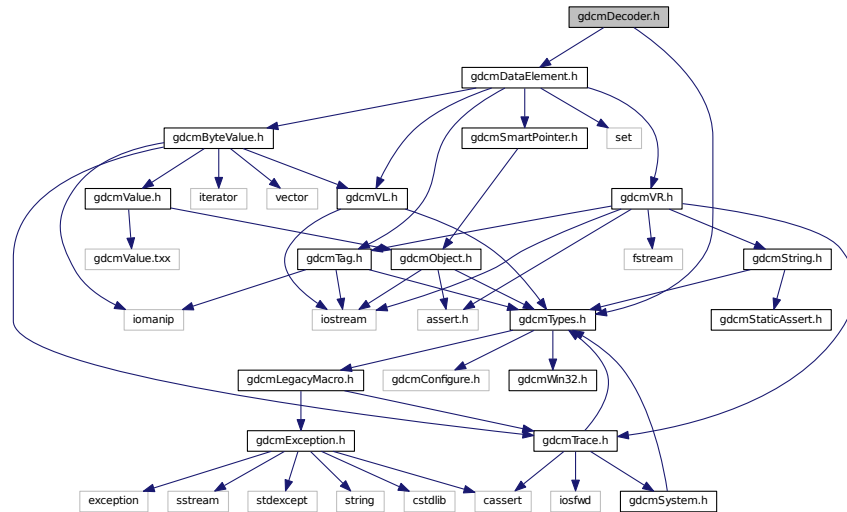

This graph shows which files directly or indirectly include this file:

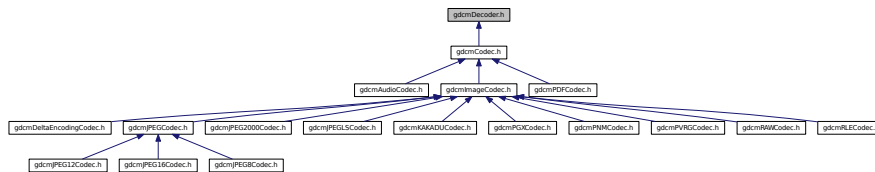

## Classes

- class `gdcm::Decoder`  
*Decoder.*

## Namespaces


- `gdcm`

## 28.61 gdcmDefinedTerms.h File Reference

```
#include "gdcmTypes.h"
```



Include dependency graph for gdcMDefinedTerms.h:

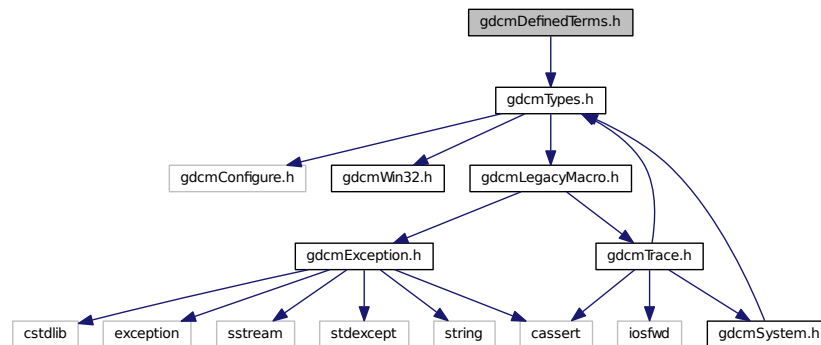

## Classes

- class [gdcM::DefinedTerms](#)

*Defined Terms are used when the specified explicit Values may be extended by implementors to include additional new Values. These new Values shall be specified in the Conformance Statement (see PS 3.2) and shall not have the same meaning as currently defined Values in this standard. A Data [Element](#) with Defined Terms that does not contain a [Value](#) equivalent to one of the Values currently specified in this standard shall not be considered to have an invalid value. Note: Interpretation [Type](#) ID (4008,0210) is an example of a Data [Element](#) having Defined Terms. It is defined to have a [Value](#) that may be one of the set of standard Values; REPORT or AMENDMENT (see PS 3.3). Because this Data [Element](#) has Defined Terms other Interpretation [Type](#) IDs may be defined by the implementor.*


## Namespaces

- [gdcM](#)

## 28.62 gdcMDeflateStream.h File Reference

```
#include "zipstreamimpl.h"
```

Include dependency graph for gdcMDeflateStream.h:

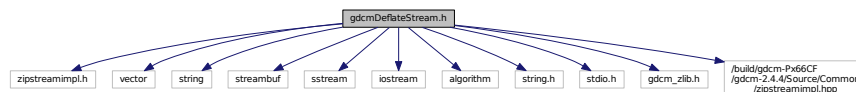

## 28.63 gdcMDefs.h File Reference


```
#include "gdcMModules.h"
```



## 28.64 gdcmDeltaEncodingCodec.h File Reference

```
#include "gdcmImageCodec.h"
```

Include dependency graph for gdcmDeltaEncodingCodec.h:

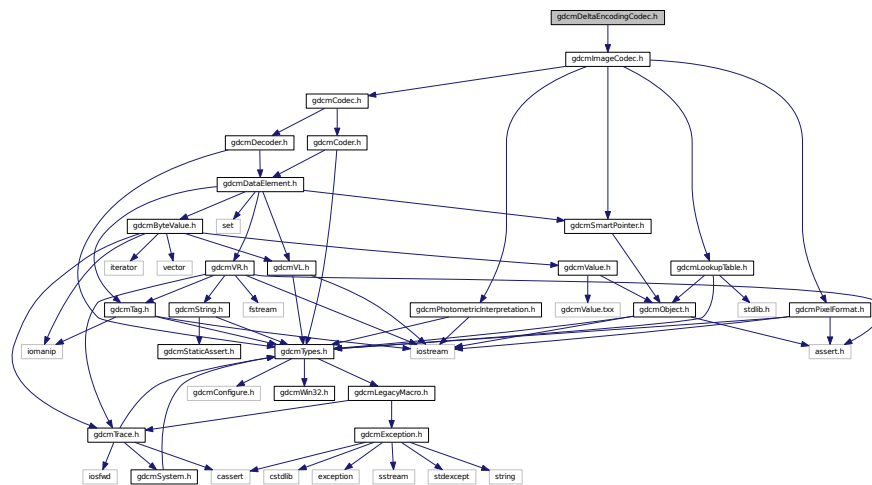

### Classes

- class [gdcm::DeltaEncodingCodec](#)

*DeltaEncodingCodec* compression used by some private vendor.


### Namespaces

- [gdcm](#)

## 28.65 gdcmDICOMDIR.h File Reference

```
#include "gdcmFileSet.h"
```

Include dependency graph for `gdcmDICOMDIR.h`:

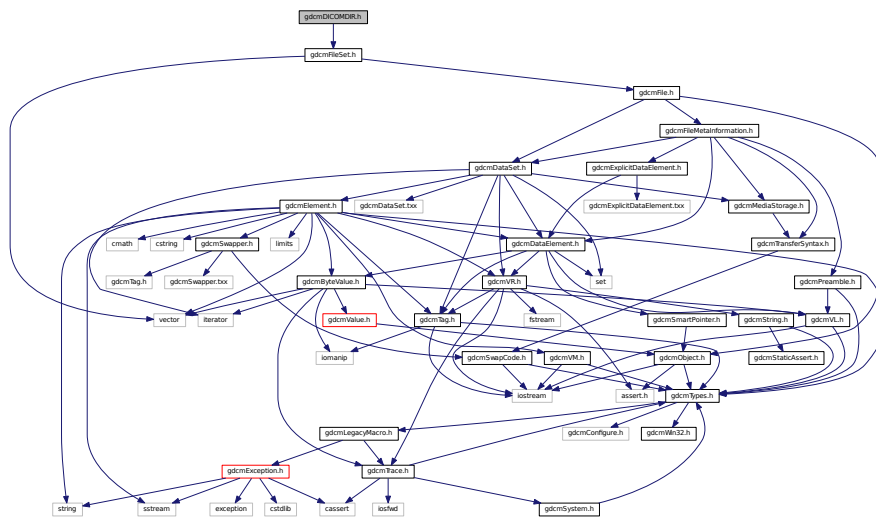

## Classes

- class `gdcm::DICOMDIR`

*DICOMDIR* class.


## Namespaces

- `gdcm`

## 28.66 gdcmDICOMDIRGenerator.h File Reference

```
#include "gdcmDirectory.h"
#include "gdcmTag.h"
#include <utility>
```

Include dependency graph for gdcmDICOMDIRGenerator.h:

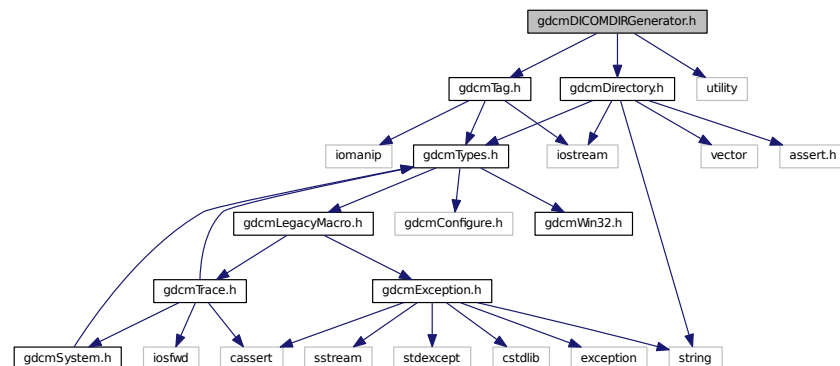

## Classes

- class [gdcm::DICOMDIRGenerator](#)

*[DICOMDIRGenerator](#) class This is a STD-GEN-CD [DICOMDIR](#) generator. ref: PS 3.11-2008 Annex D (Normative) - General Purpose CD-R and DVD Interchange Profiles.*

## Namespaces



- [gdcm](#)

## 28.67 gdcmDict.h File Reference

```

#include "gdcmTypes.h"
#include "gdcmTag.h"
#include "gdcmPrivateTag.h"
#include "gdcmDictEntry.h"
#include "gdcmSystem.h"
#include <iostream>
#include <iomanip>
#include <map>

```

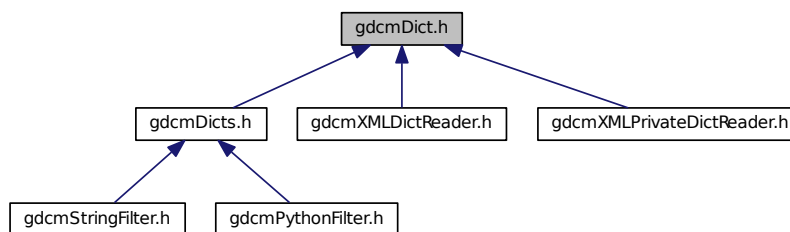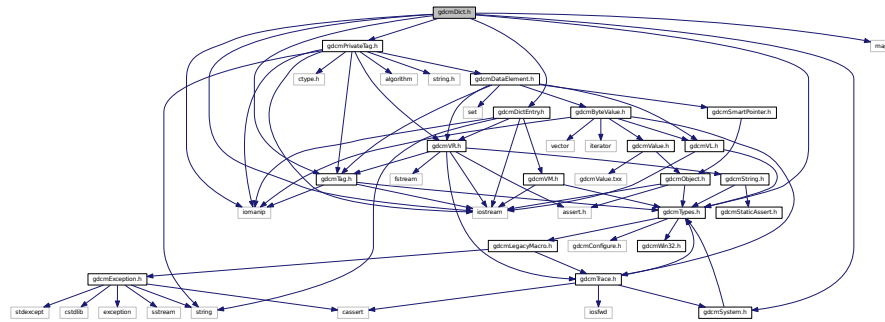

- class `gdcmm::Dict`  
*Class to represent a map of `DictEntry`.*
- class `gdcmm::PrivateDict`  
*Private `Dict`.*


- [gcdm](#)

- `std::ostream & gdcm::operator<< (std::ostream &os, const Dict &val)`
- `std::ostream & gdcm::operator<< (std::ostream &os, const PrivateDict &val)`

## 28.68 gdcmDictConverter.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmVR.h"
#include "gdcmVM.h"
#include <string>
```

Include dependency graph for gdcmDictConverter.h:

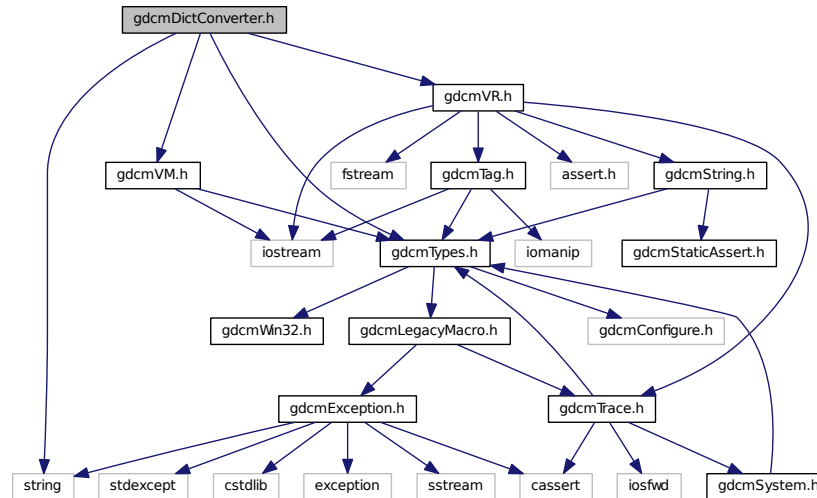

### Classes

- class [gdcm::DictConverter](#)

*Class to convert a .dic file into something else:*


### Namespaces

- [gdcm](#)


## 28.69 gdcmDictEntry.h File Reference

```
#include "gdcmVR.h"
#include "gdcmVM.h"
#include <string>
#include <iostream>
#include <iomanip>
```

Include dependency graph for `gdcmDictEntry.h`:

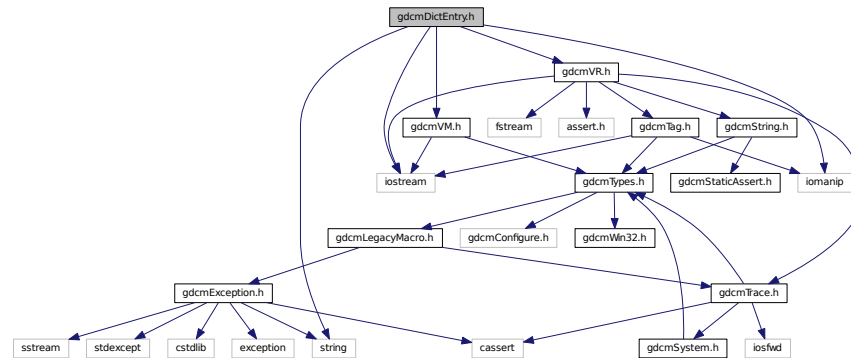

This graph shows which files directly or indirectly include this file:

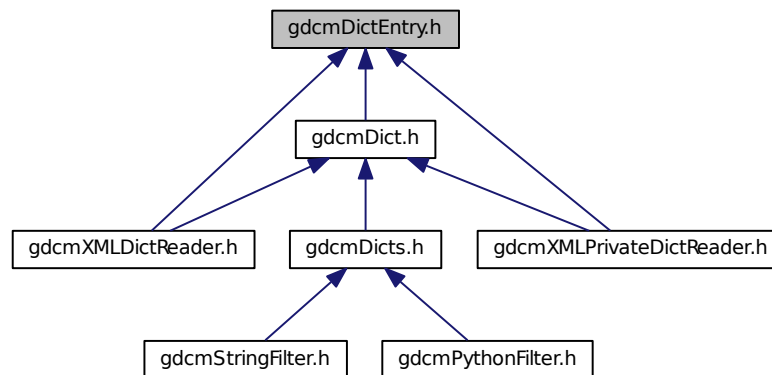

## Classes

- class [gdcm::DictEntry](#)

*Class to represent an Entry in the [Dict](#). Does not really exist within the DICOM definition, just a way to minimize storage and have a mapping from [gdcm::Tag](#) to the needed information.*

## Namespaces


- [gdcm](#)

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const DictEntry &val)`



```
#include "gdcmPrinter.h"
Include dependency graph for gdcmDictPrinter.h:
```

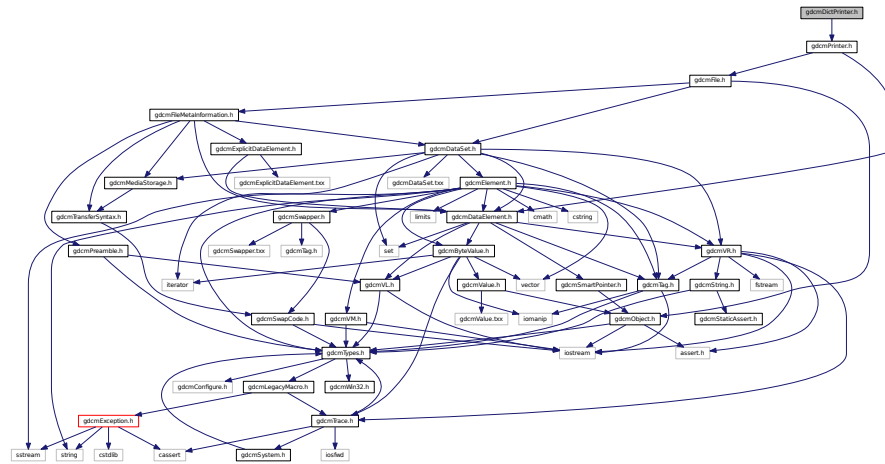


- class `gdcm::DictPrinter`

*DictPrinter* class.


- gdc

```
#include "gdcDict.h"
#include "gdcCSAHeaderDict.h"
#include <string>
```

Include dependency graph for gdcMDicts.h:

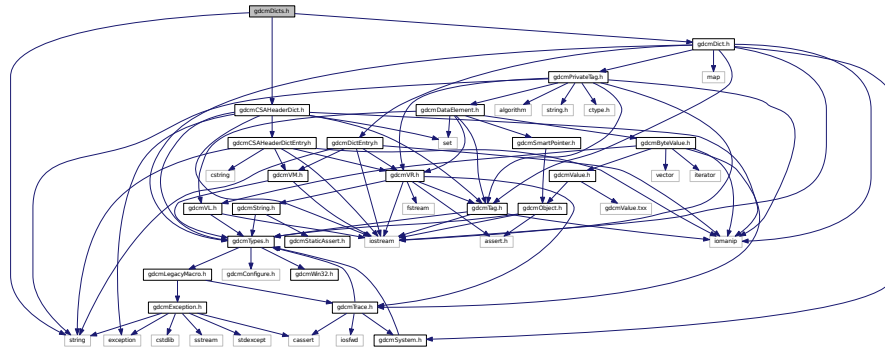

This graph shows which files directly or indirectly include this file:

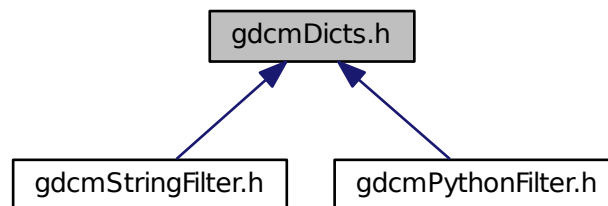

## Classes

- class `gdc::Dicts`  
*Class to manipulate the sum of knowledge (all the dict user load)*

## Namespaces

- **gdcm**

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const Dicts &d)`

## 28.72 gdcmdiff.dox File Reference

## 28.73 gdcmDIMSE.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmDIMSE.h:

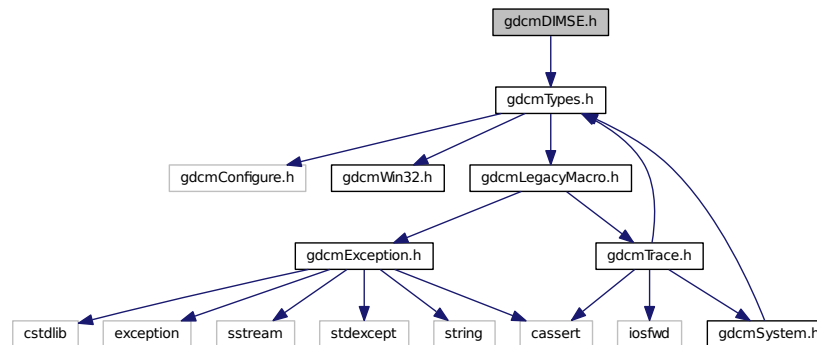

### Classes

- class [gdcm::network::CEchoRQ](#)

*[CEchoRQ](#) this file defines the messages for the cecho action.*

- class [gdcm::network::CEchoRSP](#)

*[CEchoRSP](#) this file defines the messages for the cecho action.*

- class [gdcm::network::CFind](#)

- class [gdcm::network::DIMSE](#)

*[DIMSE](#) PS 3.7 - 2009 Annex E [Command](#) Dictionary (Normative) E.1 REGISTRY OF DICOM COMMAND ELEMENTS  
[Table E.1-1](#) COMMAND FIELDS (PART 1)*


### Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.74 gdcmDirectionCosines.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmdirDirectionCosines.h`:

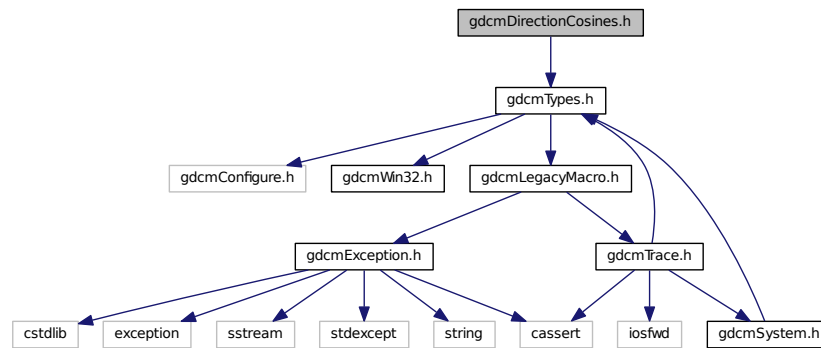

## Classes

- class `gdcmdir::DirectionCosines`  
*class to handle `DirectionCosines`*

## Namespaces

- `gdcmdir`


## 28.75 gdcmdirDirectory.h File Reference

```


#include "gdcmdirTypes.h"
#include <string>
#include <vector>
#include <iostream>
#include <assert.h>

```

Include dependency graph for `gdcmdirDirectory.h`:

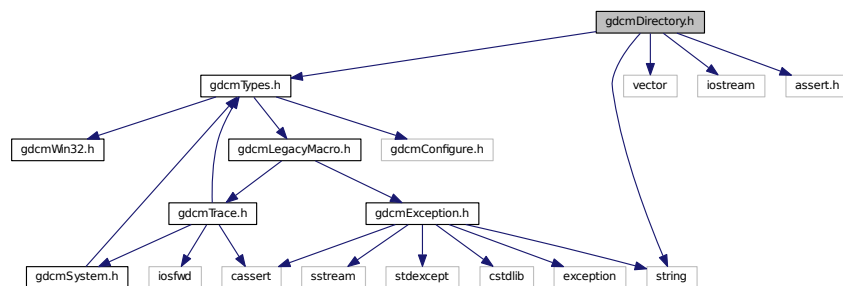

This graph shows which files directly or indirectly include this file:

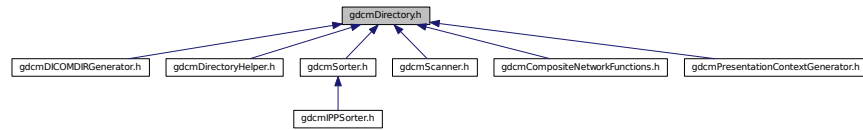

## Classes

- class [gdcm::Directory](#)  
Class for manipulation directories.

## Namespaces

- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const Directory &d)`

## 28.76 gdcmDirectoryHelper.h File Reference

```
#include "gdcmDirectory.h"
```

```
#include "gdcmDataSet.h"
```

Include dependency graph for gdcmDirectoryHelper.h:

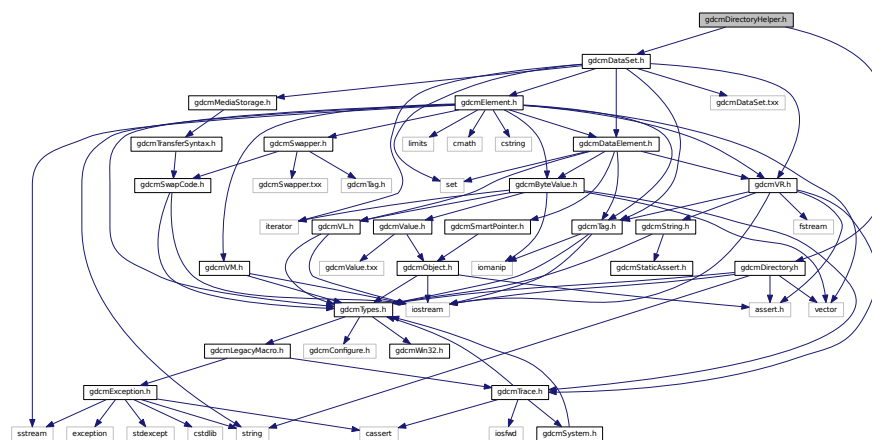

## Classes

- class [gdcm::DirectoryHelper](#)

*DirectoryHelper* this class is designed to help mitigate some of the commonly performed operations on directories. namely: 1) the ability to determine the number of series in a directory by what type of series is present 2) the ability to find all ct series in a directory 3) the ability to find all mr series in a directory 4) to load a set of DataSets from a series that's already been sorted by the IPP sorter 5) For rtstruct stuff, you need to know the sopinstanceuid of each z plane, so there's a retrieval function for that 6) then a few other functions for rtstruct writeouts.


## Namespaces

- [gdcm](#)

## 28.77 gdcmDummyValueGenerator.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmDummyValueGenerator.h:

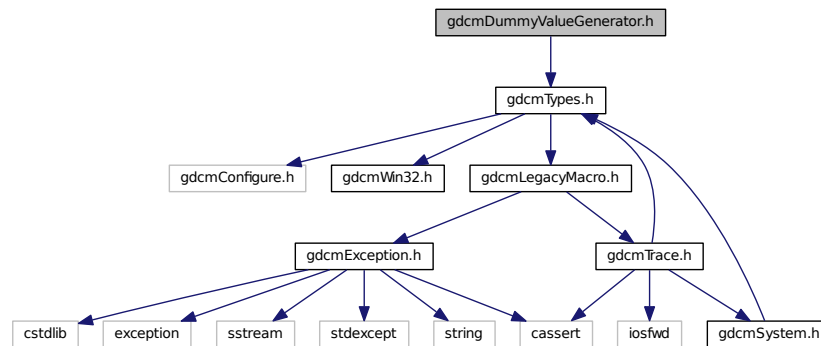

## Classes

- class [gdcm::DummyValueGenerator](#)  
Class for generating dummy value.

## Namespaces

- [gdcm](#)

## 28.78 gdcmdump.dox File Reference

## 28.79 gdcmDumper.h File Reference

```
#include "gdcmPrinter.h"
```

- class `gdcm::Dumper`

## Namespaces

- ## 28.80 gdcmElement.h File Reference

Generated on Sat Sep 26 2015 07:03:31 for GDCM by Doxygen






## Functions

- ignore\_char const [gdcm::backslash](#) ("\\")
- std::istream & [gdcm::operator>>](#) (std::istream &in, ignore\_char const &ic)
- template<typename Float >  
std::string [gdcm::to\\_string](#) (Float data)

## 28.81 gdcmEncapsulatedDocument.h File Reference

```
#include "gdcmFile.h"
```

Include dependency graph for gdcmEncapsulatedDocument.h:

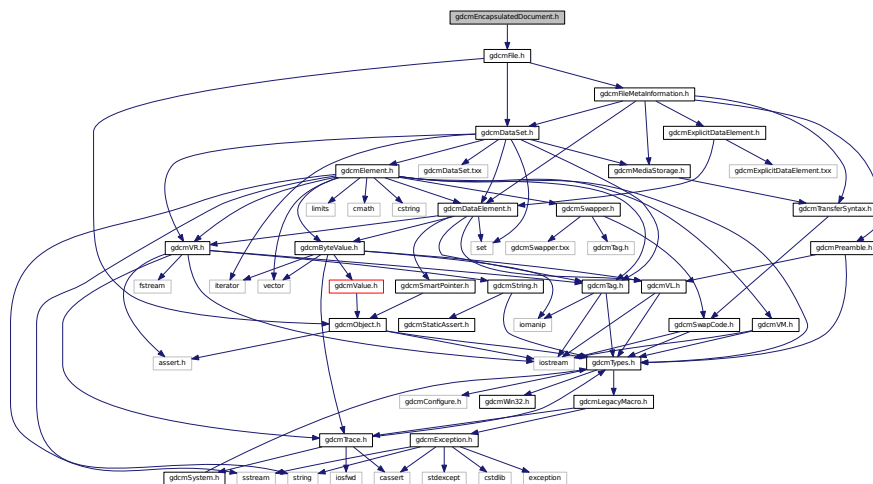

## Classes

- class [gdcm::EncapsulatedDocument](#)  
*EncapsulatedDocument.*


## Namespaces

- [gdcm](#)

## 28.82 gdcmEnumeratedValues.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmEnumeratedValues.h`:

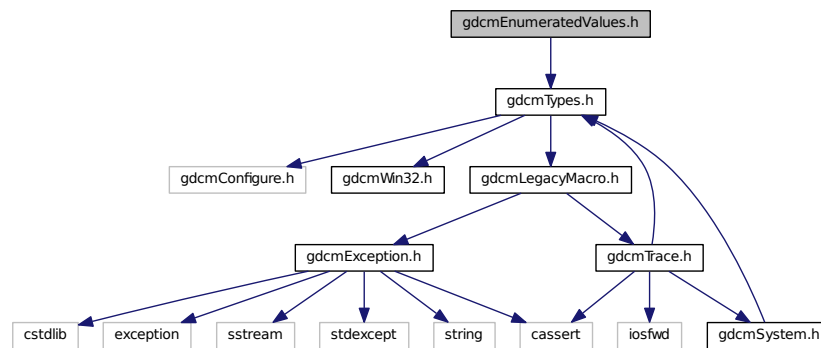

## Classes

- class [gdcm::EnumeratedValues](#)

***Element.** A Data **Element** with Enumerated Values that does not have a **Value** equivalent to one of the Values specified in this standard has an invalid value within the scope of a specific Information Object/SOP Class definition. Note:*


## Namespaces

- [gdcm](#)


## 28.83 gdcmEvent.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmEvent.h`:

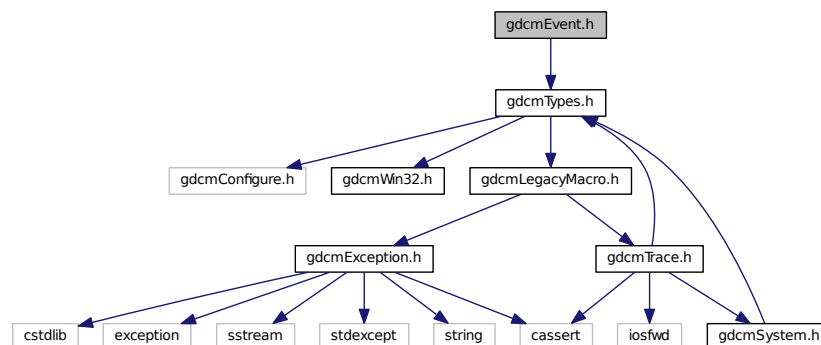

This graph shows which files directly or indirectly include this file:

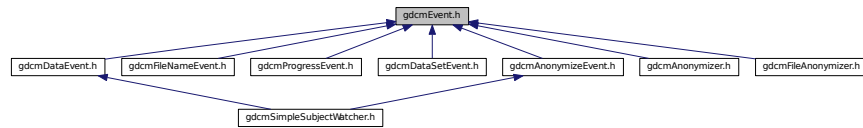

## Classes

- class [gdcm::AbortEvent](#)
- class [gdcm::AnyEvent](#)
- class [gdcm::EndEvent](#)
- class [gdcm::Event](#)  
*superclass for callback/observer methods*
- class [gdcm::ExitEvent](#)
- class [gdcm::InitializeEvent](#)
- class [gdcm::IterationEvent](#)
- class [gdcm::ModifiedEvent](#)
- class [gdcm::NoEvent](#)
- class [gdcm::StartEvent](#)
- class [gdcm::UserEvent](#)

## Namespaces

- [gdcm](#)

## Macros

- `#define gdcmEventMacro(classname, super)`

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, Event &e)`  
*Generic inserter operator for [Event](#) and its subclasses.*

### 28.83.1 Macro Definition Documentation

#### 28.83.1.1 `#define gdcmEventMacro( classname, super )`

**Value:**

```

\
class classname : public super { \
    public: \
        typedef classname Self; \
        typedef super Superclass; \
        classname() {} \
        virtual ~classname() {} \

```

```

virtual const char * GetEventName() const { return #classname; } \
virtual bool CheckEvent(const ::gdcM::Event* e) const \
{ return dynamic_cast<const Self*>(e) ? true : false; } \
virtual ::gdcM::Event* MakeObject() const \
{ return new Self; } \
classname(const Self&s) : super(s){}; \
private: \
void operator=(const Self&); \
}

```


## 28.84 gdcMException.h File Reference

```

#include <cassert>
#include <cstdlib>
#include <exception>
#include <sstream>
#include <stdexcept>
#include <string>

```

Include dependency graph for gdcMException.h:

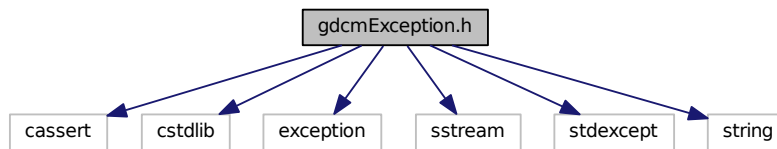

This graph shows which files directly or indirectly include this file:



### Classes

- class [gdcM::Exception](#)  
*Exception.*

### Namespaces

- [gdcM](#)


## 28.85 gdcMExplicitDataElement.h File Reference

```


#include "gdcMDataElement.h"
#include "gdcMExplicitDataElement.txx"

```

Include dependency graph for gdcmExplicitDataElement.h:

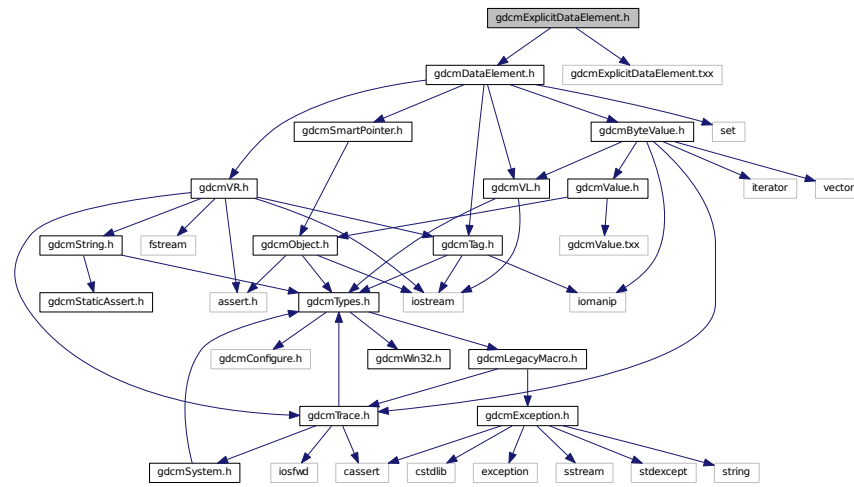

This graph shows which files directly or indirectly include this file:

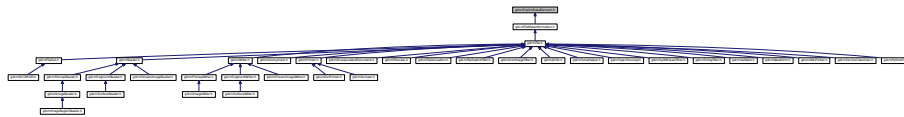

## Classes

- class [gdcm::ExplicitDataElement](#)

*Class to read/write a [DataElement](#) as Explicit Data [Element](#).*


## Namespaces

- [gdcm](#)

## 28.86 gdcmExplicitImplicitDataElement.h File Reference

```
#include "gdcmDataElement.h"
#include "gdcmExplicitImplicitDataElement.txx"
```

Include dependency graph for `gdcmExplicitImplicitDataElement.h`:

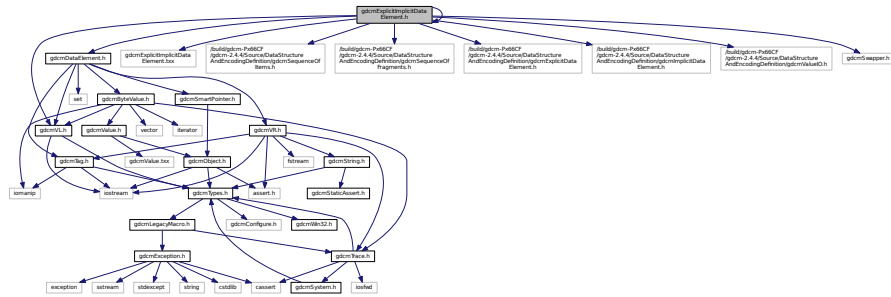

## Classes

- class `gdcm::ExplicitImplicitDataElement`  
*Class to read/write a `DataElement` as ExplicitImplicit Data `Element`.*


## Namespaces

- **gdcm**

## 28.87 gdcmFiducials.h File Reference

```
#include "gdcmFile.h"
```

Include dependency graph for gdcMfiducials.h:

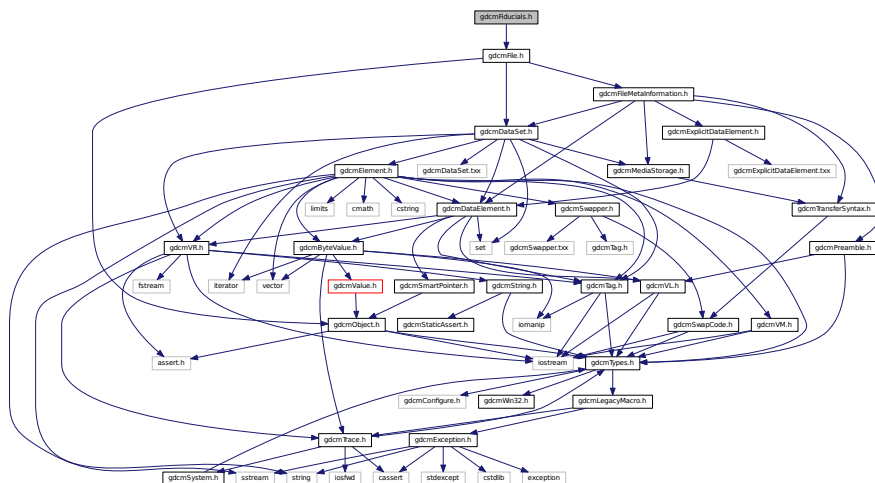


## Classes

- class `gdcm::Fiducials`


## Namespaces

- gdcm

```
#include "gdcmObject.h"
#include "gdcmDataSet.h"
#include "gdcmFileMetaInformation.h"
Include dependency graph for gdcmFile.h:
```

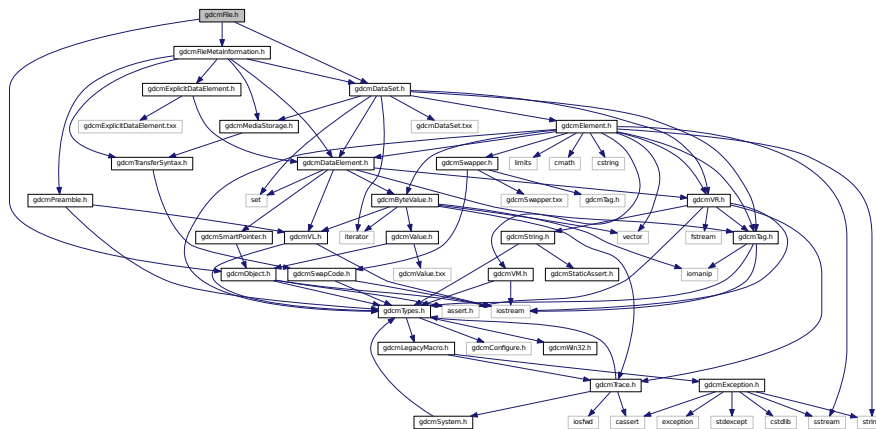

This graph shows which files directly or indirectly include this file:

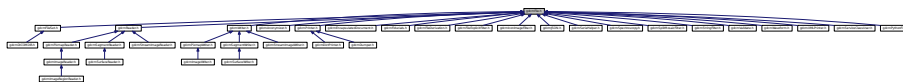

## Classes

- class `gdcm::File`

a DICOM [File](#) See PS 3.10 [File](#): A [File](#) is an ordered string of zero or more bytes, where the first byte is at the beginning of the file and the last byte at the end of the [File](#). Files are identified by a unique [File](#) ID and may be written, read and/or deleted.

## Namespaces

- `gdcm`


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const File &val)`

## 28.89 gdcmFileAnonymizer.h File Reference

```
#include "gdcmSubject.h"
#include "gdcmEvent.h"
#include "gdcmTag.h"
#include "gdcmVL.h"
```

Include dependency graph for `gdcmFileAnonymizer.h`:

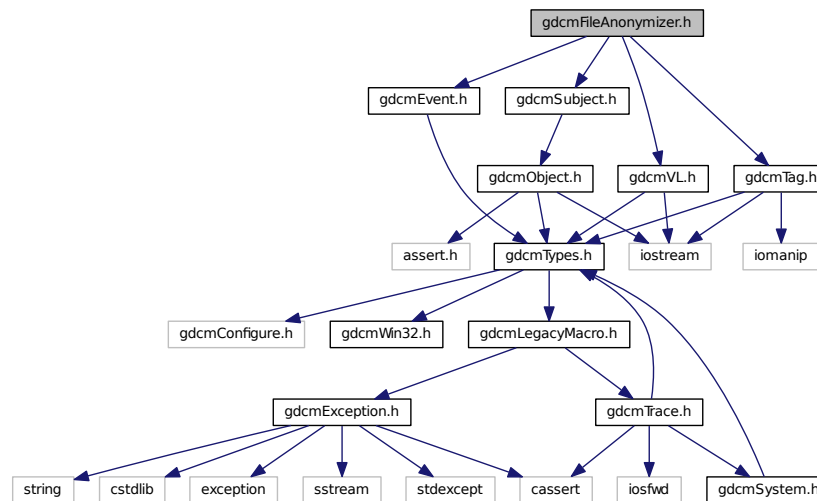

## Classes

- class `gdcm::FileAnonymizer`  
*[FileAnonymizer](#)*.

## Namespaces


- `gdcm`

## 28.90 gdcmFileChangeTransferSyntax.h File Reference

```
#include "gdcmSubject.h"
#include "gdcmSmartPointer.h"
```



Include dependency graph for gdcmFileChangeTransferSyntax.h:

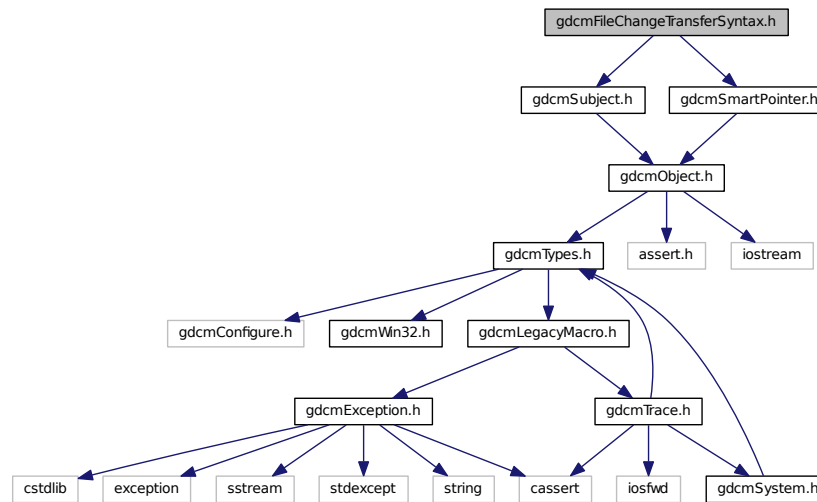

## Classes

- class [gdcm::FileChangeTransferSyntax](#)

*FileChangeTransferSyntax.*


## Namespaces

- [gdcm](#)

## 28.91 gdcmFileDerivation.h File Reference

```
#include "gdcmFile.h"
```

Include dependency graph for `gdcmFileDerivation.h`:

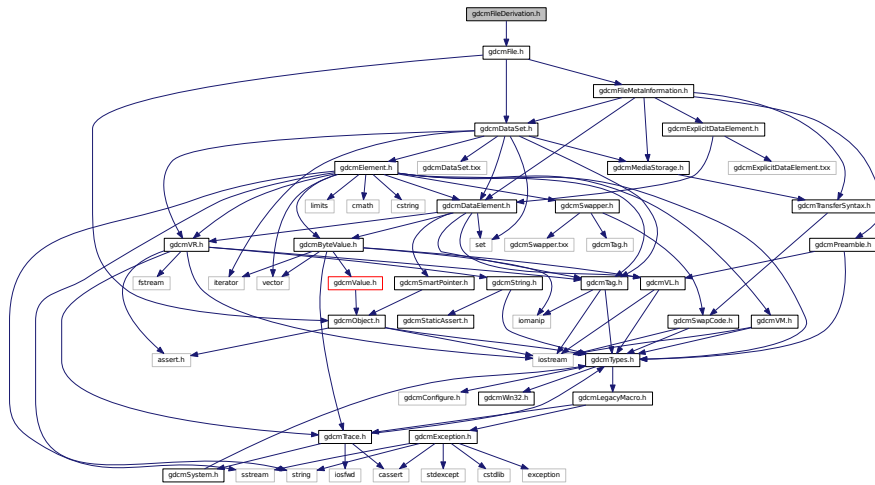

## Classes

- class `gdcm::FileDerivation`

*FileDerivation* class See PS 3.16 - 2008 For the list of Code Value that can be used for in Derivation Code Sequence.

## Namespaces


- `gdcm`

## 28.92 gdcmFileExplicitFilter.h File Reference


```
#include "gdcmFile.h"
```



Include dependency graph for `gdcmlFileMetaInformation.h`:

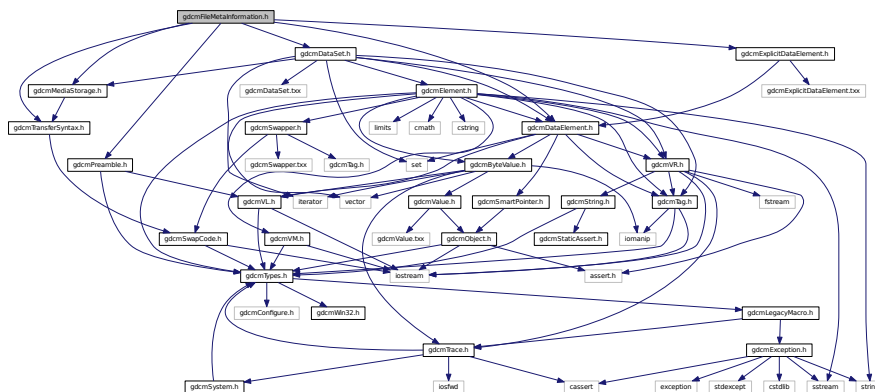

This graph shows which files directly or indirectly include this file:

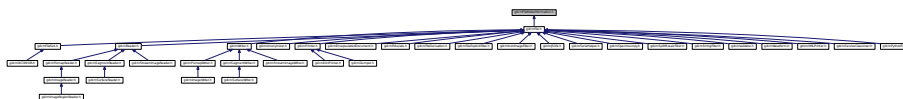

## Classes

- class [gdcml::FileMetaInformation](#)  
Class to represent a [File](#) Meta Information.

## Namespaces

- [gdcml](#)


## Functions

- `std::ostream & gdcml::operator<< (std::ostream &os, const FileMetaInformation &val)`

## 28.94 gdcmlFilename.h File Reference

```
#include "gdcmlTypes.h"
#include <string>
```

Include dependency graph for gdcmFilename.h:

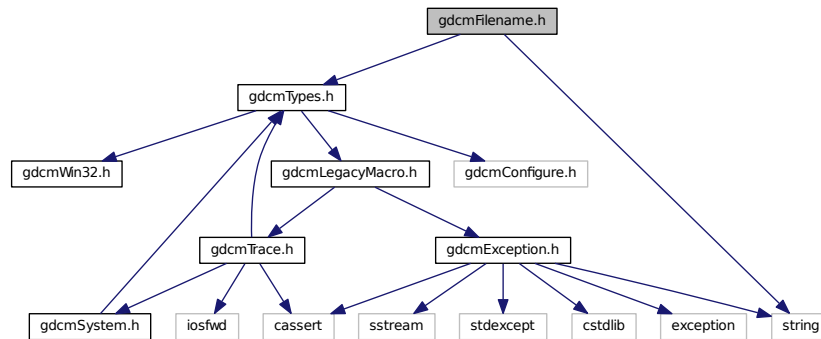

## Classes

- class [gdcm::Filename](#)  
Class to manipulate file name's.

## Namespaces


- [gdcm](#)

## 28.95 gdcmFileNameEvent.h File Reference

```
#include "gdcmEvent.h"
```

```
#include "gdcmTag.h"
```

Include dependency graph for gdcmFileNameEvent.h:

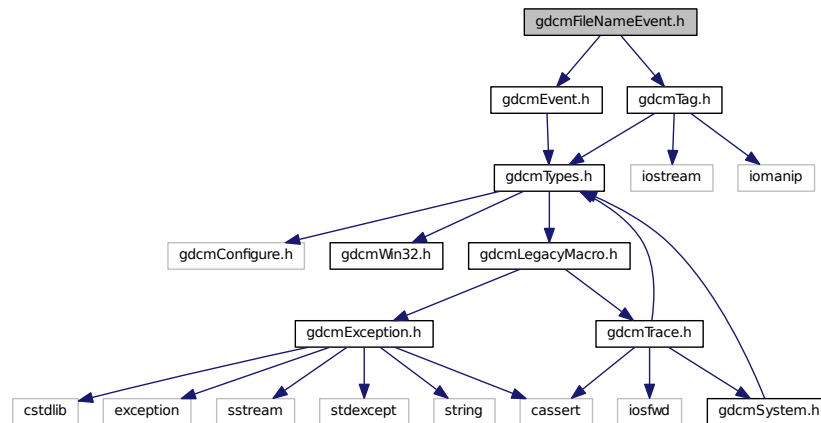

## Classes

- class [gdcm::FileNameEvent](#)  
*FileNameEvent* Special type of event triggered during processing of *FileSet*.


## Namespaces

- [gdcm](#)

## 28.96 gdcmFilenameGenerator.h File Reference

```
#include "gdcmTypes.h"
#include <string>
#include <vector>
```

Include dependency graph for gdcmFilenameGenerator.h:

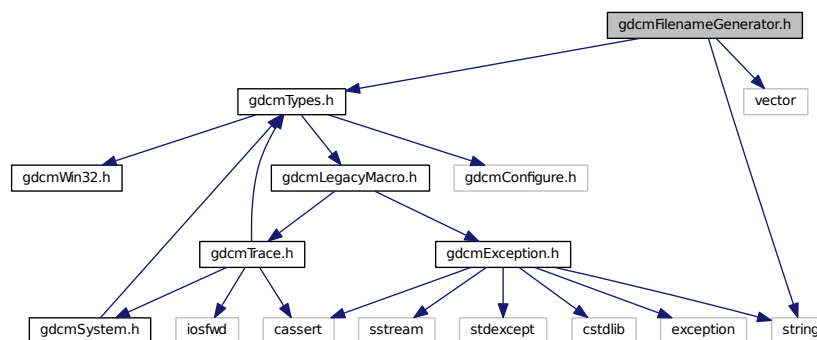

## Classes

- class [gdcm::FilenameGenerator](#)  
*FilenameGenerator*.

## Namespaces

- [gdcm](#)

## 28.97 gdcmFileSet.h File Reference

```
#include "gdcmFile.h"
#include <vector>
```

[illegible]

```
graph BT
    gdcmDICOMDIR.h --> gdcmFileSet.h
```

- class `gdcm::FileSet`


## Namespaces

- ## Functions

- Generated on Sat Sep 26 2015 07:03:31 for GDCM by Doxygen

## 28.98 gdcmFileStreamer.h File Reference

```
#include "gdcmSubject.h"
#include "gdcmSmartPointer.h"
Include dependency graph for gdcmFileStreamer.h:
```

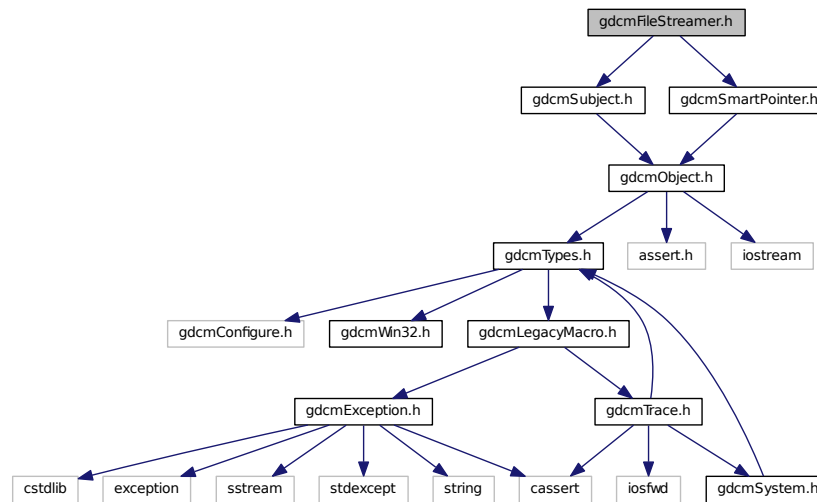

## Classes

- class [gdcm::FileStreamer](#)

*[FileStreamer](#) This class let a user create a massive DICOM [DataSet](#) from a template DICOM file, by appending chunks of data.*

## Namespaces


- [gdcm](#)

## 28.99 gdcmFindPatientRootQuery.h File Reference


```
#include "gdcmBaseRootQuery.h"
```



Include dependency graph for gdcMFindPatientRootQuery.h:

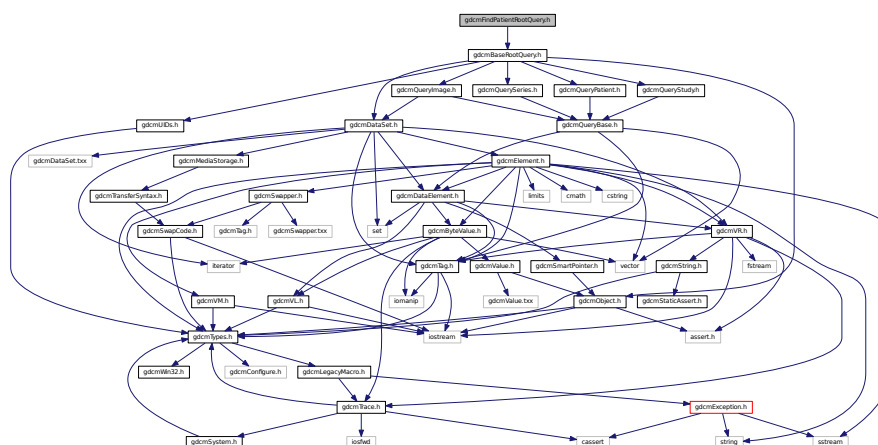

This graph shows which files directly or indirectly include this file:

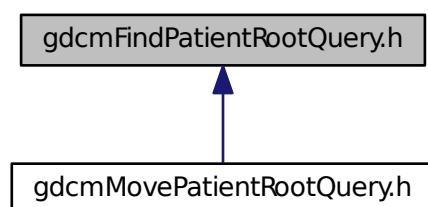

## Classes

- class `gdcm::FindPatientRootQuery`


*PatientRootQuery* contains: the class which will produce a dataset for c-find with patient root.

## Namespaces

- **gdcm**


## 28.100 gdcMFindStudyRootQuery.h File Reference

```
#include "gdcmBaseRootQuery.h"
```

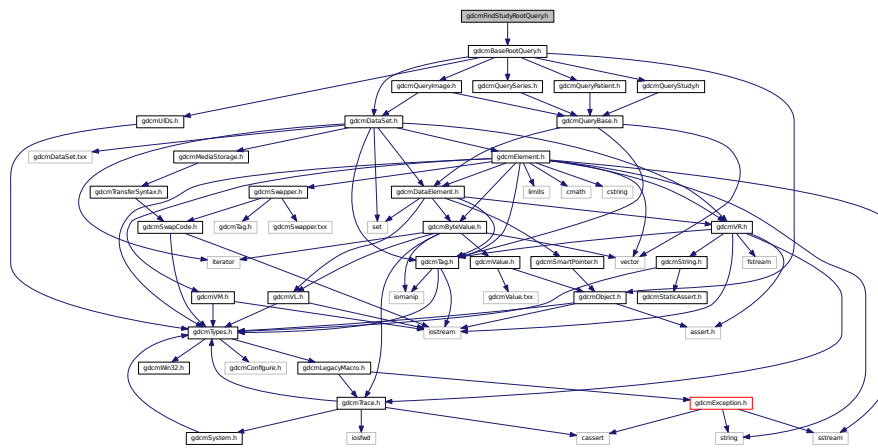


```
#include "gdcmDataElement.h"
#include "gdcmByteValue.h"
#include "gdcmSmartPointer.h"
#include "gdcmParseException.h"
```

Include dependency graph for gdcmFragment.h:

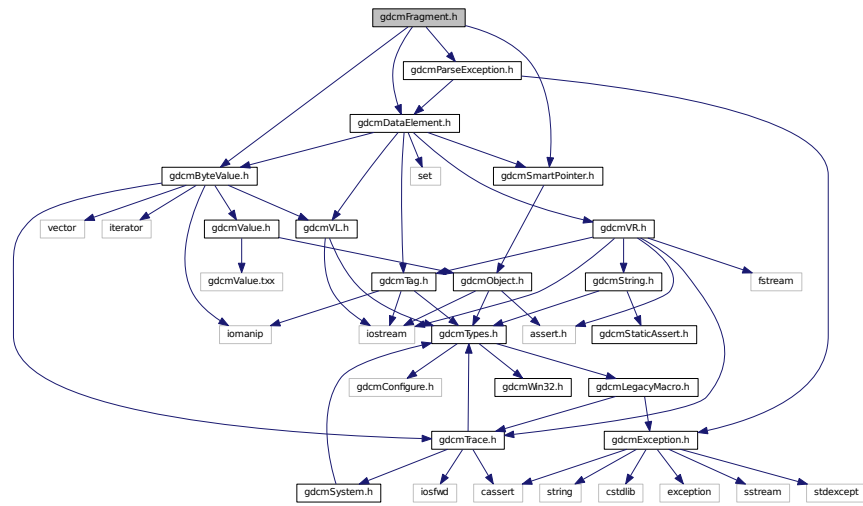

This graph shows which files directly or indirectly include this file:

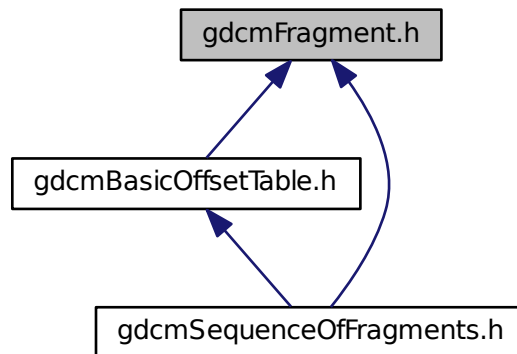

## Classes

- class [gdcm::Fragment](#)  
Class to represent a *Fragment*.

## Namespaces

- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const Fragment &val)`

## 28.102 [gdcmgendir.dox](#) File Reference

## 28.103 [gdcmGlobal.h](#) File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmGlobal.h`:

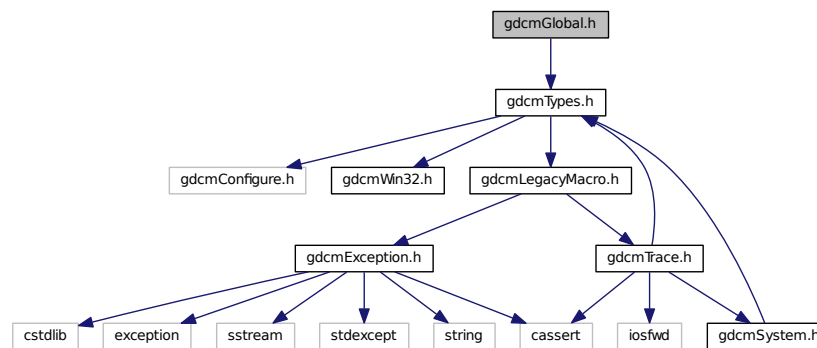

## Classes

- class `gdcm::Global`  
*Global.*

## Namespaces

- `gdcm`

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Global &g)`


## Variables

- static Global `gdcm::GlobalInstance`

## 28.104 gdcmGroupDict.h File Reference

```
#include "gdcmTypes.h"
#include <assert.h>
#include <vector>
#include <string>
#include <iostream>
#include <iomanip>
```

Include dependency graph for gdcmGroupDict.h:

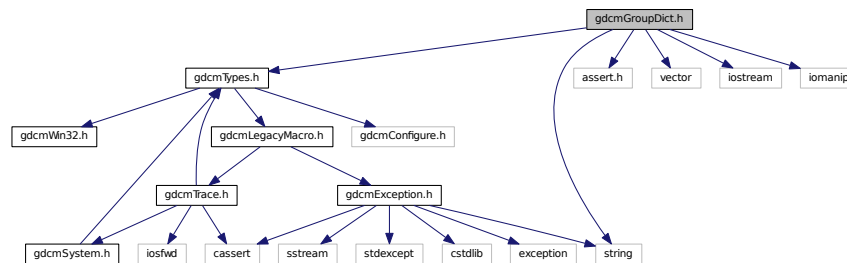

### Classes

- class [gdcm::GroupDict](#)

*Class to represent the mapping from group number to its abbreviation and name.*

### Namespaces

- [gdcm](#)


### Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const GroupDict &_val)`


## 28.105 gdcmIconImage.h File Reference

```
#include "gdcmBitmap.h"
```

Include dependency graph for `gdcmIconImage.h`:

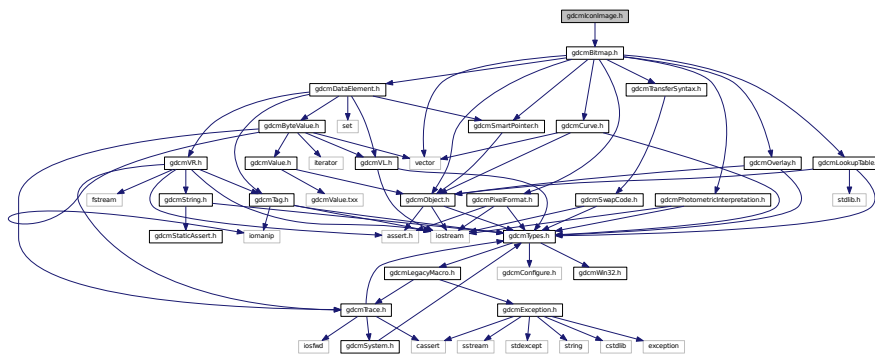

This graph shows which files directly or indirectly include this file:

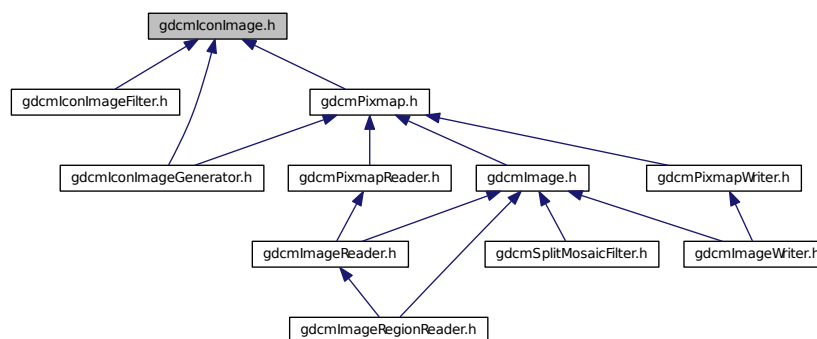

## Namespaces

- `gdcm`

## Typedefs

- typedef Bitmap `gdcm::IconImage`

## 28.106 gdcmIconImageFilter.h File Reference

```
#include "gdcmFile.h"
#include "gdcmIconImage.h"
```


[illegible]

- class `gdcm::IconImageFilter`

## Namespaces

- ## 28.107 gdcmlconImageGenerator.h File Reference

Include dependency graph for `gdcmIconImageGenerator.h`:

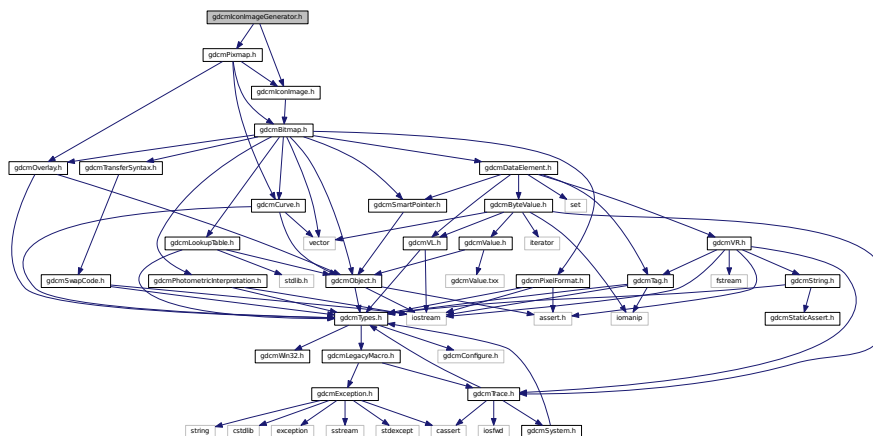





## Classes

- class [gdcm::Image](#)

*Image* This is the container for an *Image* in the general sense. From this container you should be able to request information like:


## Namespaces

- [gdcm](#)

## 28.109 gdcmImageApplyLookupTable.h File Reference

```
#include "gdcmImageToImageFilter.h"
```

Include dependency graph for gdcmImageApplyLookupTable.h:

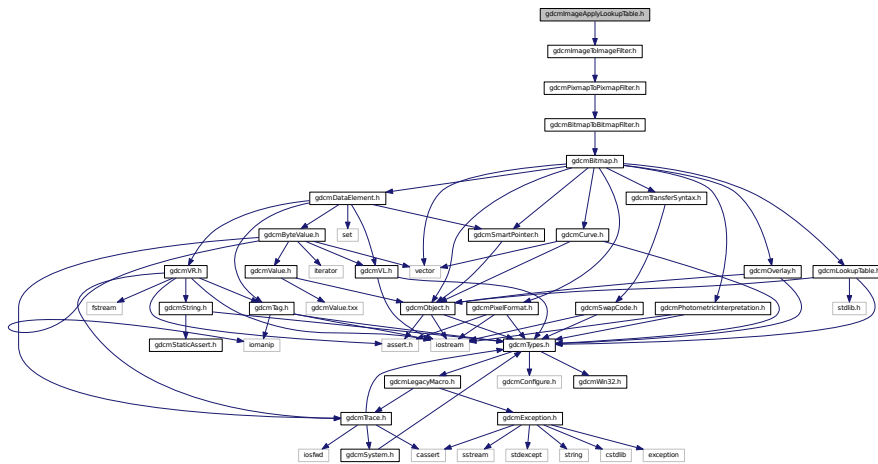

## Classes

- class [gdcm::ImageApplyLookupTable](#)

*ImageApplyLookupTable* class It applies the LUT the *PixelData* (only *PALETTE\_COLOR* images) Output will be a *PhotometricInterpretation=RGB* image.

## Namespaces


- [gdcm](#)

## 28.110 gdcmImageChangePhotometricInterpretation.h File Reference

```
#include "gdcmImageToImageFilter.h"
```

```
#include "gdcmPhotometricInterpretation.h"
```

Include dependency graph for `gdcmImageChangePhotometricInterpretation.h`:

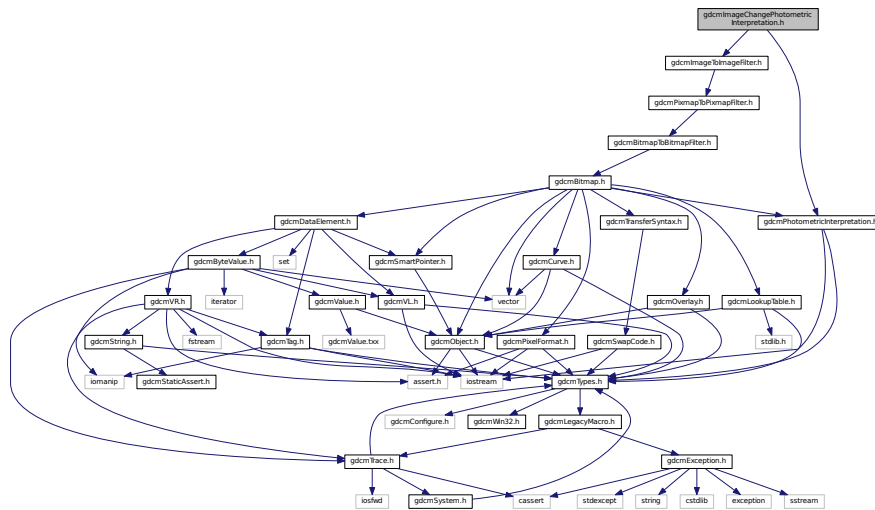

## Classes

- class [gdcm::ImageChangePhotometricInterpretation](#)

*ImageChangePhotometricInterpretation* class Class to change the Photometric Interpretation of an input DICOM.


## Namespaces

- [gdcm](#)

## 28.111 gdcmImageChangePlanarConfiguration.h File Reference

```
#include "gdcmImageToImageFilter.h"
```

Include dependency graph for gdcmImageChangePlanarConfiguration.h:

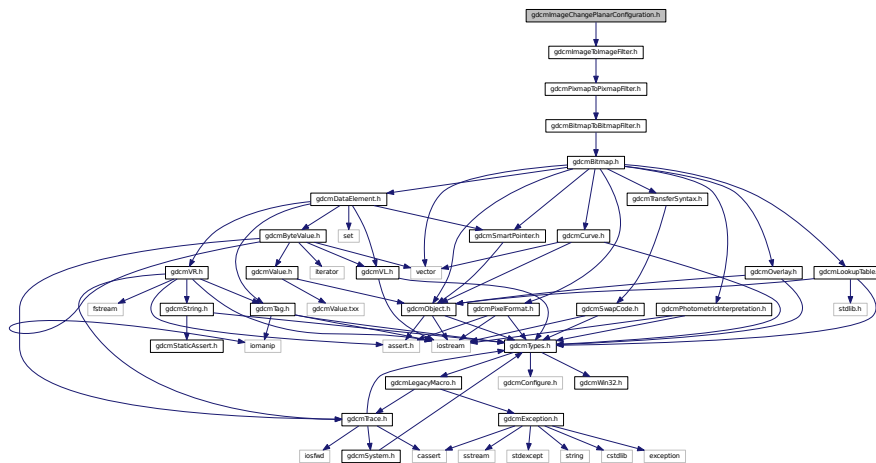

## Classes

- class [gdcm::ImageChangePlanarConfiguration](#)

*[ImageChangePlanarConfiguration](#) class Class to change the Planar configuration of an input DICOM By default it will change into the more usual representation: PlanarConfiguration = 0.*

## Namespaces

- [gdcm](#)

## 28.112 gdcmImageChangeTransferSyntax.h File Reference

```
#include "gdcmImageToImageFilter.h"
#include "gdcmTransferSyntax.h"
```




- ImageChangeTransferSyntax** class Class to change the transfer syntax of an input DICOM.

- **gdcm**

```
#include "gdcmCodec.h"
#include "gdcmPhotometricInterpretation.h"
#include "gdcmLookupTable.h"
#include "gdcmSmartPointer.h"
#include "gdcmPixelFormat.h"
```

Include dependency graph for gdcmImageCodec.h:

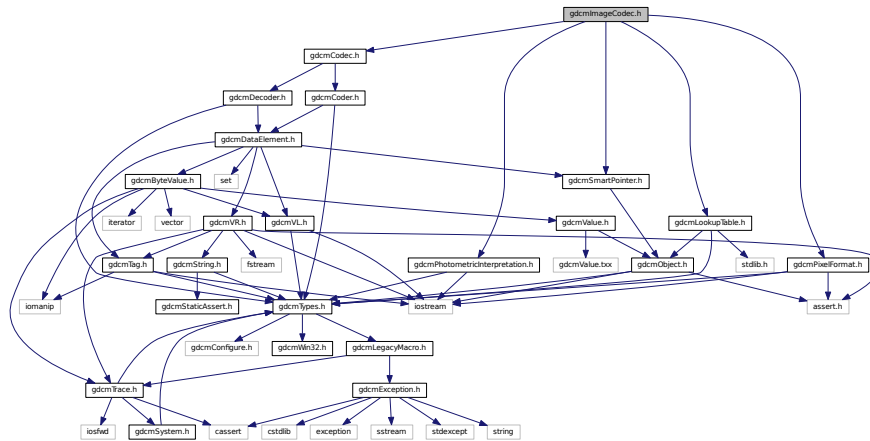

This graph shows which files directly or indirectly include this file:



## Classes

- class [gdcm::ImageCodec](#)


*ImageCodec.*

## Namespaces

- [gdcm](#)

## 28.114 gdcmImageConverter.h File Reference


```
#include "gdcmTypes.h"
```

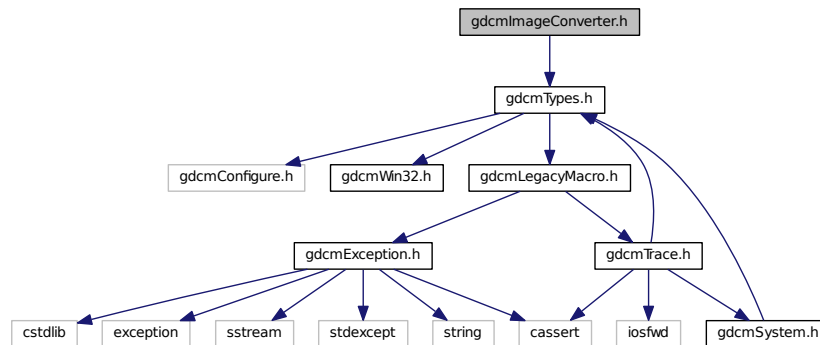

- class `gdcm::ImageConverter`  
*Image Converter.*

- **gdc**

Include dependency graph for `gdcmlImageFragmentSplitter.h`:

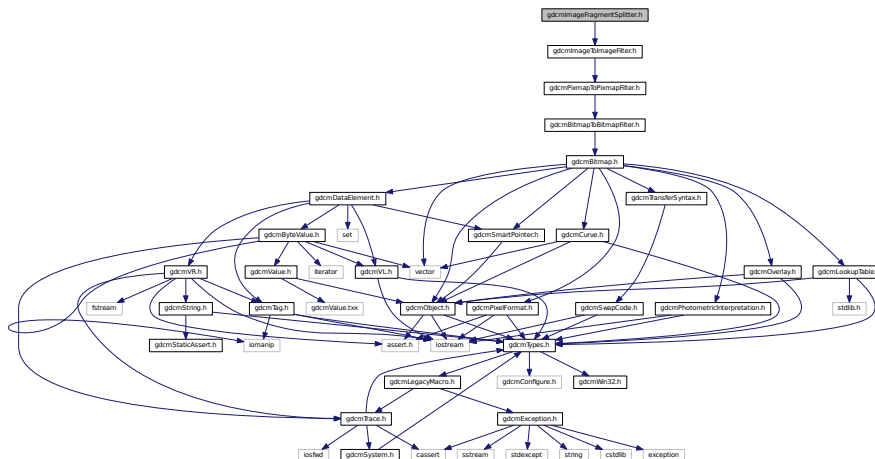

## Classes

- class [gdcm::ImageFragmentSplitter](#)

*ImageFragmentSplitter* class For single frame image, DICOM standard allow splitting the frame into multiple fragments.


## Namespaces

- [gdcm](#)

## 28.116 gdcmImageHelper.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmTag.h"
#include <vector>
#include "gdcmPixelFormat.h"
#include "gdcmPhotometricInterpretation.h"
#include "gdcmSmartPointer.h"
#include "gdcmLookupTable.h"
```

Include dependency graph for gdcmImageHelper.h:

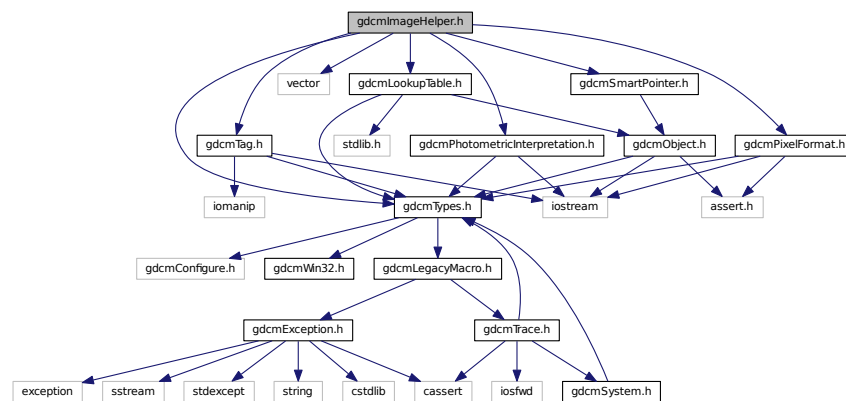

## Classes

- class [gdcm::ImageHelper](#)

*ImageHelper* (internal class, not intended for user level)


## Namespaces

- [gdcm](#)





```
#include "gdcmImage.h"
#include "gdcmRegion.h"
```

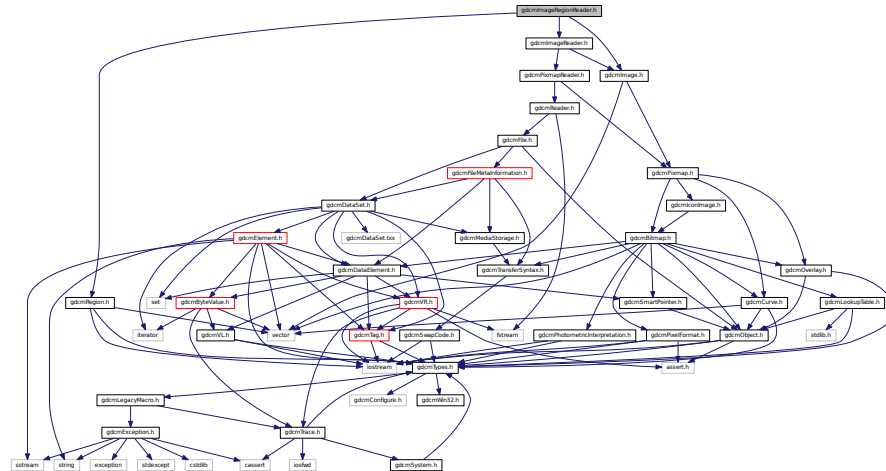

## Classes

- class `gdcm::ImageRegionReader`

*ImageRegionReader.*

## Namespaces

- **gdcm**

## 28.119 gdcmlImageToImageFilter.h File Reference

```
#include "gdcmPixmapToPixmapFilter.h"
```




- class `gdcm::ImageWriter`  
*ImageWriter*.


- `gdcm`

## 28.122 gdcmImplementationClassUIDSub.h File Reference

Include dependency graph for `gdcmlImplementationClassUIDSub.h`:

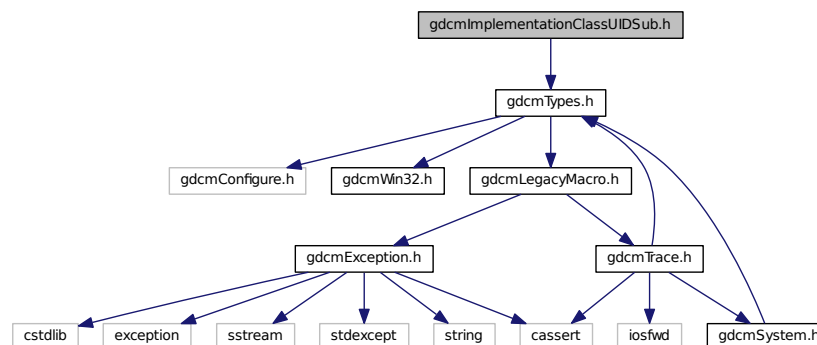

This graph shows which files directly or indirectly include this file:

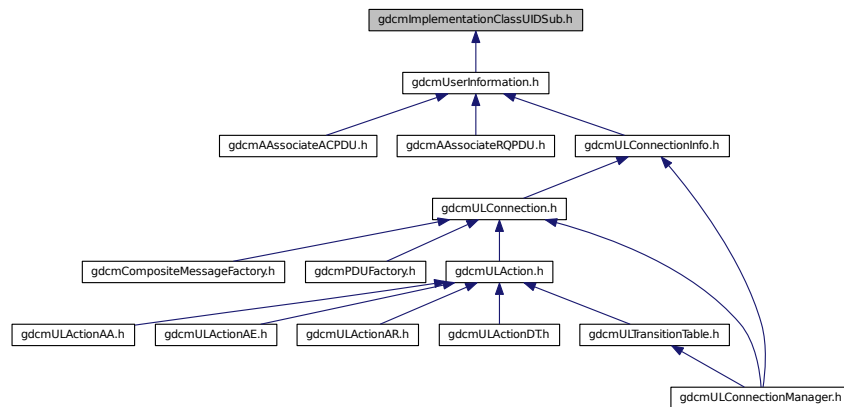

## Classes

- class [gdcml::network::ImplementationClassUIDSub](#)

[ImplementationClassUIDSub](#) PS 3.7 [Table D.3-1](#) IMPLEMENTATION CLASS UID SUB-ITEM FIELDS (A-ASSOCIAT↔E-RQ)


## Namespaces

- [gdcml](#)
- [gdcml::network](#)

## 28.123 gdcmlImplementationUIDSub.h File Reference

```
#include "gdcmlTypes.h"
```

Include dependency graph for gdcmlImplementationUIDSub.h:

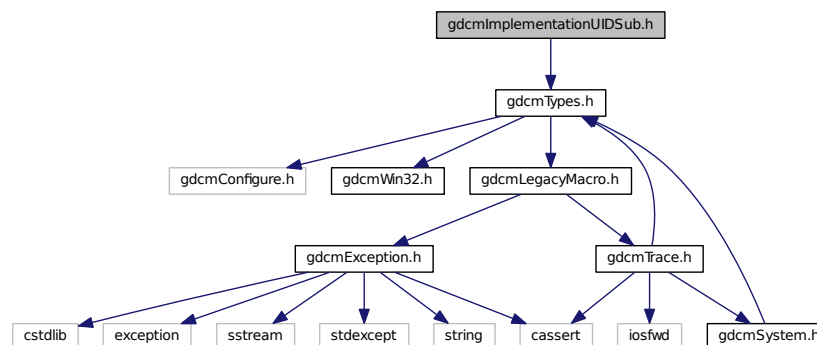

## Classes

- class [gdcm::network::ImplementationUIDSub](#)

*ImplementationUIDSub Table D.3-2 IMPLEMENTATION UID SUB-ITEM FIELDS (A-ASSOCIATE-AC)*


## Namespaces

- [gdcm](#)
- [gdcm::network](#)


## 28.124 gdcmImplementationVersionNameSub.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmImplementationVersionNameSub.h:

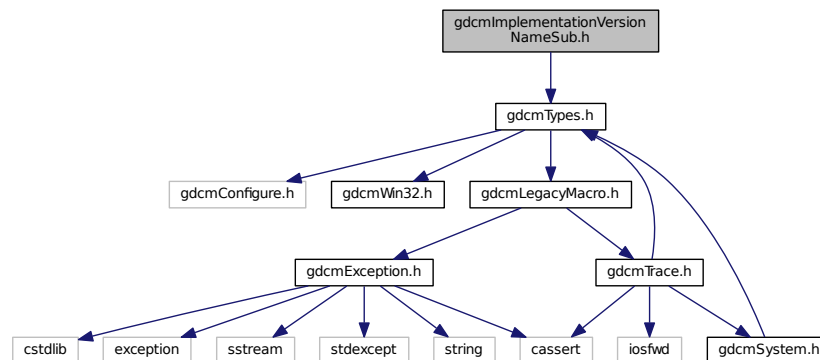

This graph shows which files directly or indirectly include this file:

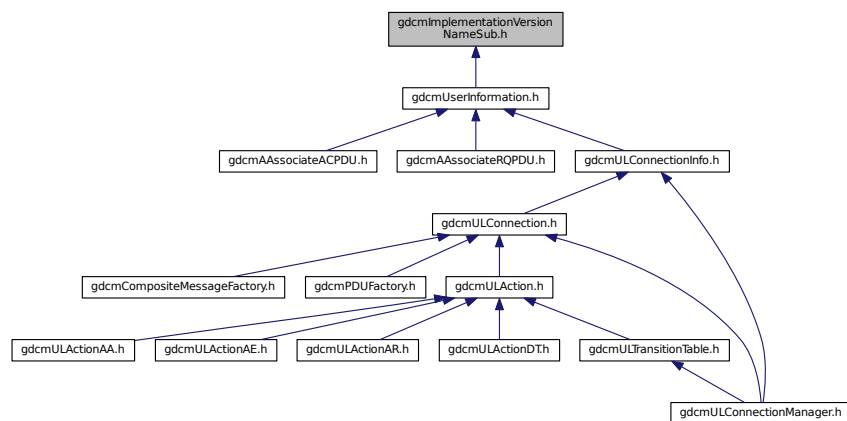

## Classes

- class [gdcm::network::ImplementationVersionNameSub](#)

*ImplementationVersionNameSub Table D.3-3 IMPLEMENTATION VERSION NAME SUB-ITEM FIELDS (A-ASSOCIATE-RQ)*

## Namespaces


- [gdcm](#)
- [gdcm::network](#)

## 28.125 gdcmImplicitDataElement.h File Reference

```
#include "gdcmDataElement.h"
```

```
#include "gdcmImplicitDataElement.txx"
```

Include dependency graph for `gdcmImplicitDataElement.h`:

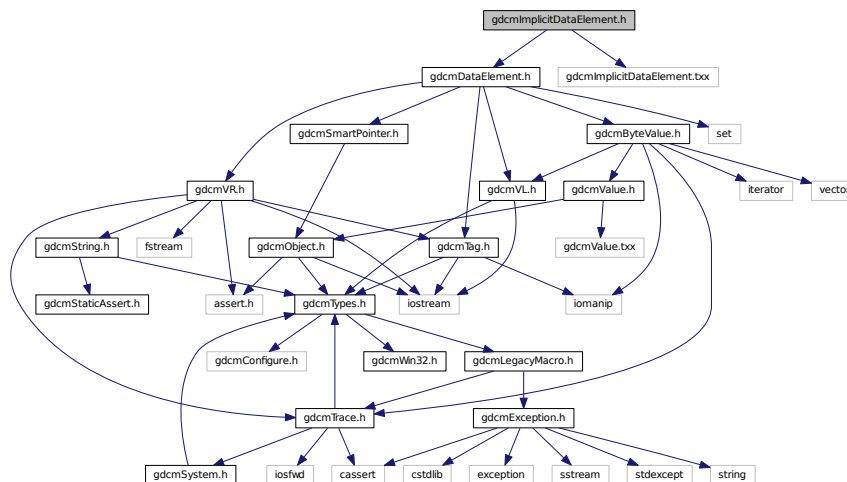

## Classes

- class [gdcm::ImplicitDataElement](#)

*Class to represent an Implicit VR Data Element.*


## Namespaces

- [gdcm](#)


## 28.126 gdcminfo.dox File Reference

## 28.127 gdcmIOD.h File Reference

```
#include "gdcmTypes.h"  
#include "gdcmTag.h"  
#include "gdcmIODEntry.h"  
#include <vector>  
Include dependency graph for gdcmIOD.h:
```

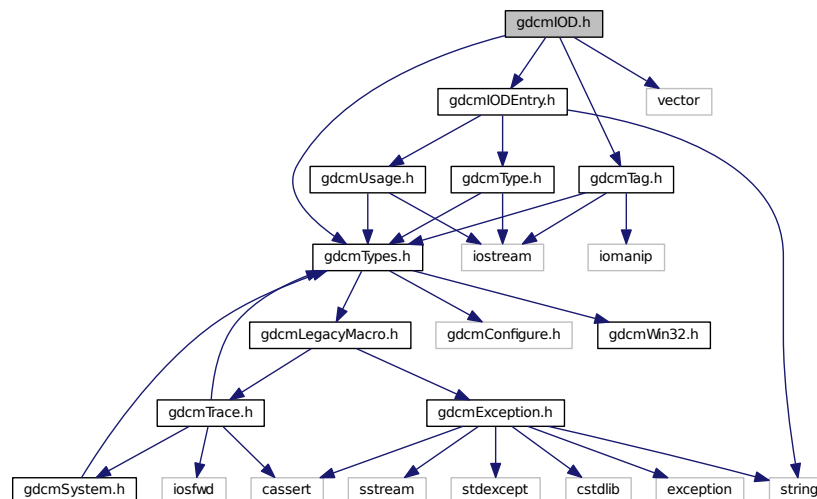

This graph shows which files directly or indirectly include this file:

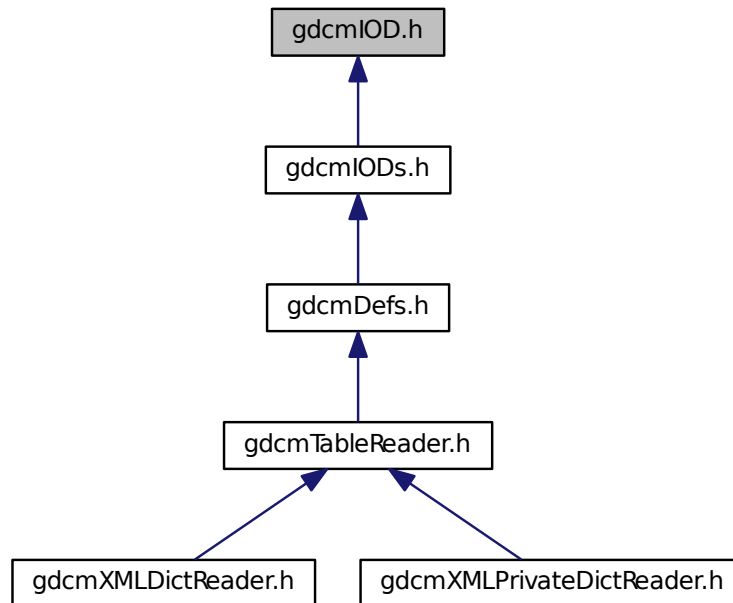

## Classes

- class [gdcml::IOD](#)  
*Class for representing a [IOD](#).*

## Namespaces

- [gdcml](#)

## Functions


- `std::ostream & gdcml::operator<< (std::ostream &_os, const IOD &_val)`

## 28.128 gdcmlIOEntry.h File Reference


```
#include "gdcmlUsage.h"
#include "gdcmlType.h"
#include <string>
```



Include dependency graph for gdcmIODEntry.h:

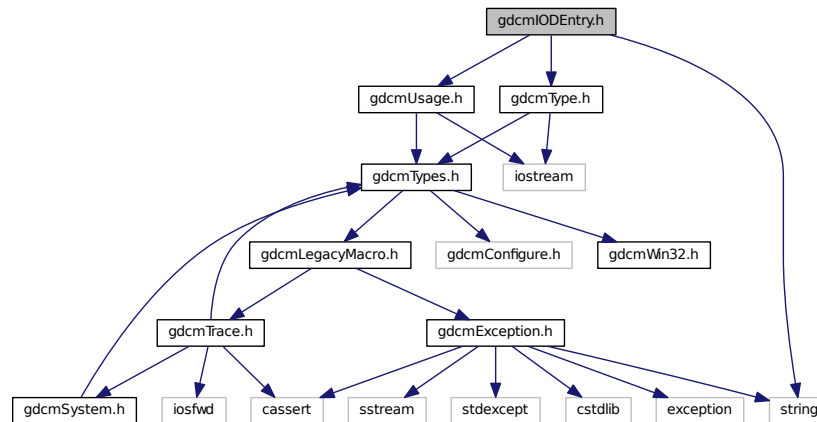

This graph shows which files directly or indirectly include this file:

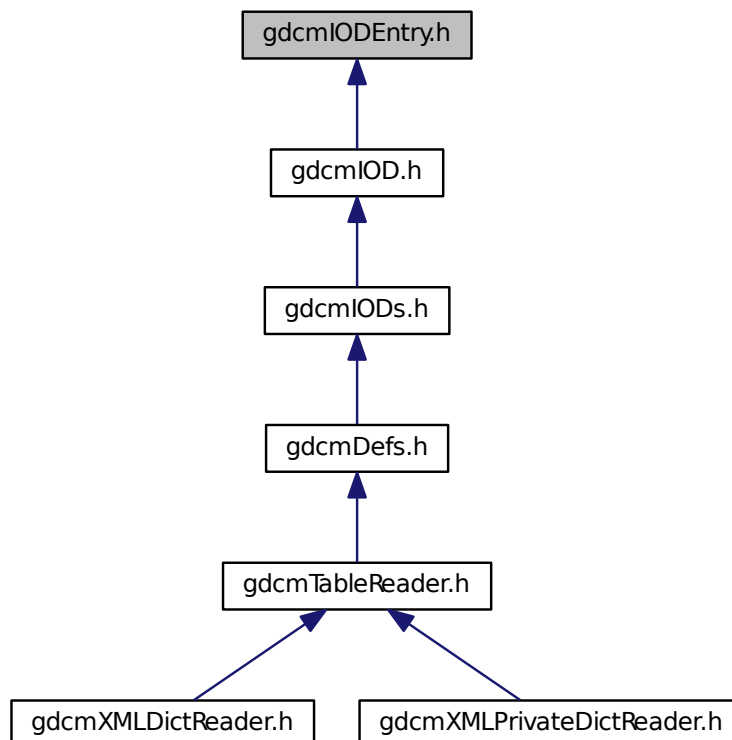

## Classes

- class [gdcm::IODEntry](#)

*Class for representing a [IODEntry](#).*

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const IODEntry &_val)`


## 28.129 gdcmIODs.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmIOD.h"
#include <map>
```

Include dependency graph for gdcmIODs.h:

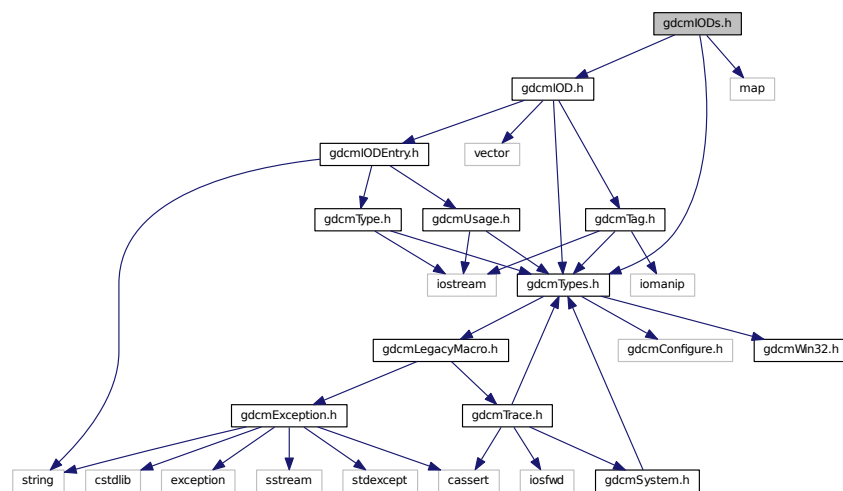

This graph shows which files directly or indirectly include this file:

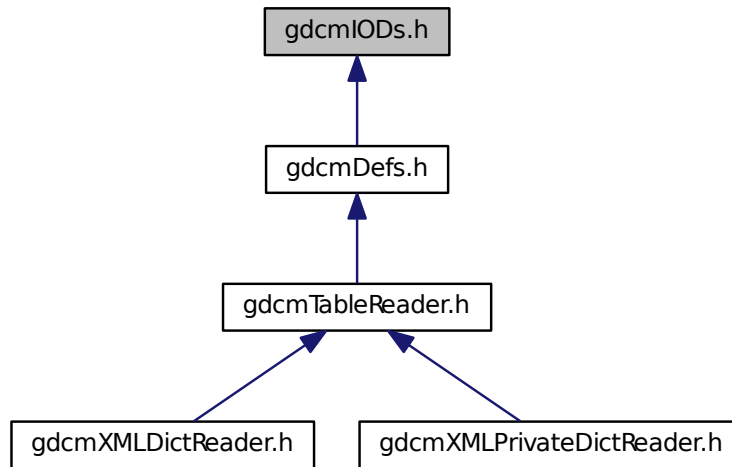

## Classes

- class [gdcm::IODs](#)

*Class for representing a [IODs](#).*

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const IODs &_val)`

## 28.130 gdcmIPPSorter.h File Reference

```
#include "gdcmSorter.h"  
#include <vector>  
#include <string>
```

Include dependency graph for `gdcmIPPSorter.h`:

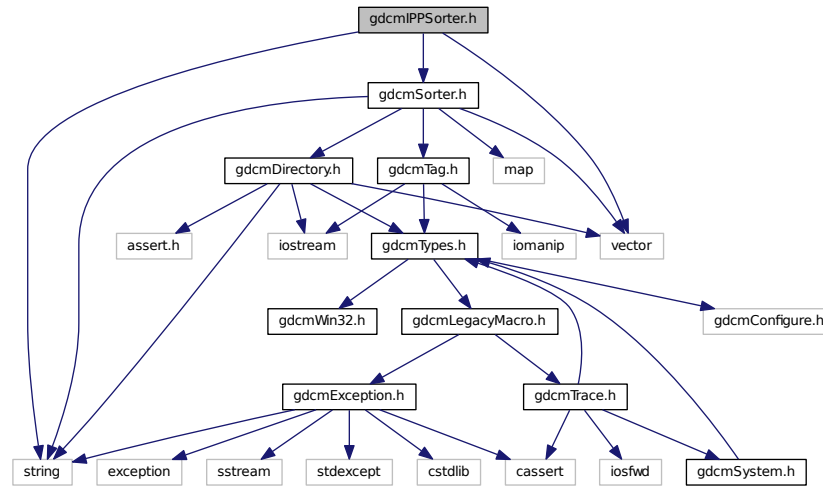

## Classes

- class `gdcm::IPPSorter`

*IPPSorter* Implement a simple *Image* Position (*Patient*) sorter, along the *Image Orientation* (*Patient*) direction. This algorithm does NOT support duplicate and will FAIL in case of duplicate IPP.

## Namespaces

- `gdcm`

## 28.131 gdcmItem.h File Reference

```

#include "gdcmDataElement.h"
#include "gdcmDataSet.h"
#include "gdcmParseException.h"
#include "gdcmSwapper.h"
#include "gdcmItem.txx"

```

[illegible]

```
graph BT; gdcmlItem.h --> gdcmlSequenceOfItems.h
```

- class `gdc::Item`

## Namespaces


- ## Functions

- Generated on Sat Sep 26 2015 07:03:31 for GDCM by Doxygen

## 28.132 gdcMJPEG12Codec.h File Reference

```
#include "gdcmJPEGCodec.h"
```

Include dependency graph for gdcMJPEG12Codec.h:

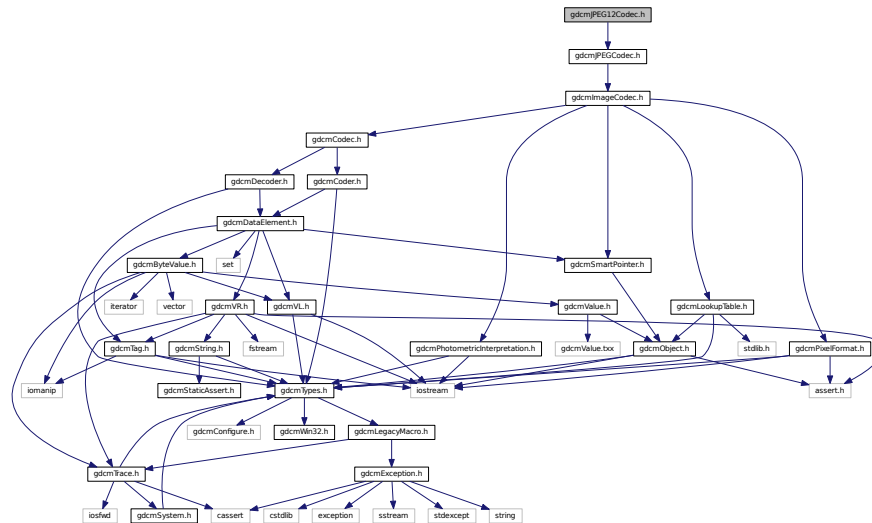

## Classes

- class `gdcm::JPEG12Codec`

*Class to do JPEG 12bits (lossy & lossless)*


## Namespaces

- **gdcm**

## 28.133 gdcMJPEG16Codec.h File Reference

```
#include "gdcmJPEGCodec.h"
```

Include dependency graph for gdcmJPEG16Codec.h:

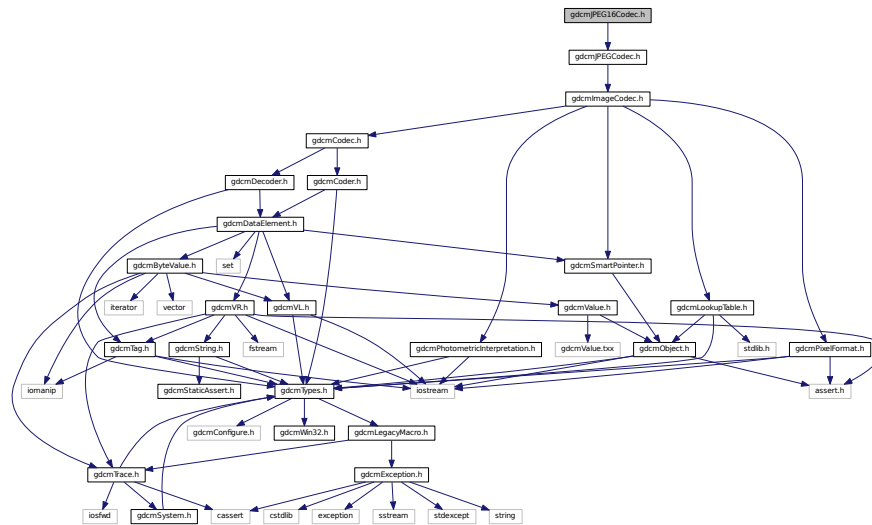

## Classes

- class [gdcm::JPEG16Codec](#)

*Class to do JPEG 16bits (lossless)*


## Namespaces

- [gdcm](#)

## 28.134 gdcmJPEG2000Codec.h File Reference

```
#include "gdcmImageCodec.h"
```

Include dependency graph for `gdcMJPEG2000Codec.h`:

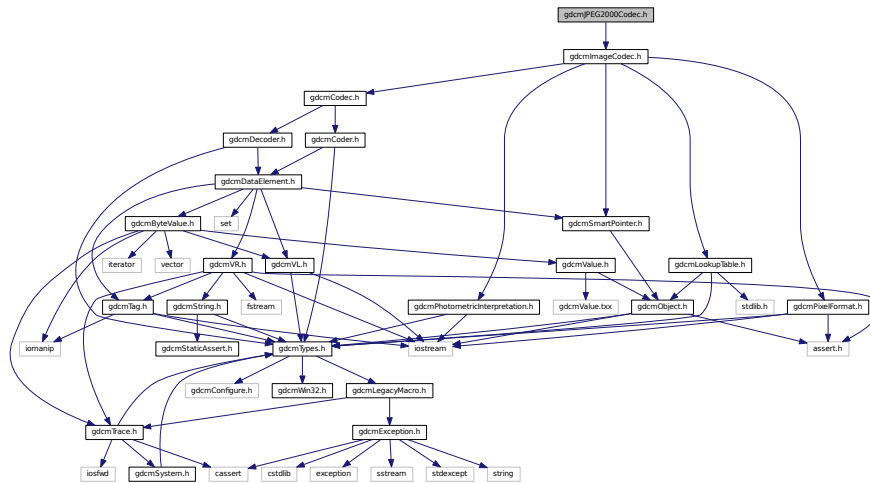

## Classes

- class [gdcM::JPEG2000Codec](#)

*Class to do JPEG 2000.*

## Namespaces


- [gdcM](#)

## 28.135 gdcMJPEG8Codec.h File Reference

```
#include "gdcMJPEGCodec.h"
```



Include dependency graph for gdcmJPEG8Codec.h:

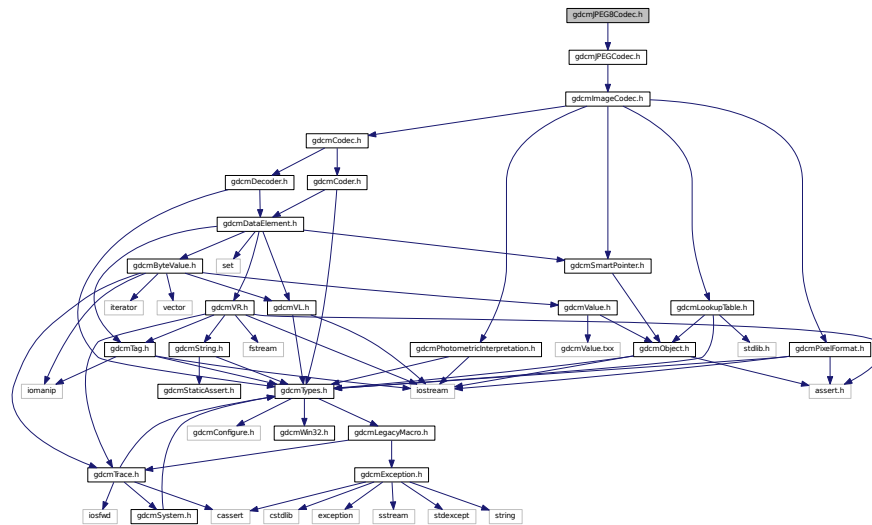

## Classes

- class [gdcm::JPEG8Codec](#)

*Class to do JPEG 8bits (lossy & lossless)*


## Namespaces

- [gdcm](#)


## 28.136 gdcmJPEGCodec.h File Reference

```
#include "gdcmImageCodec.h"
```

Include dependency graph for `gdcmJPEGCodec.h`:

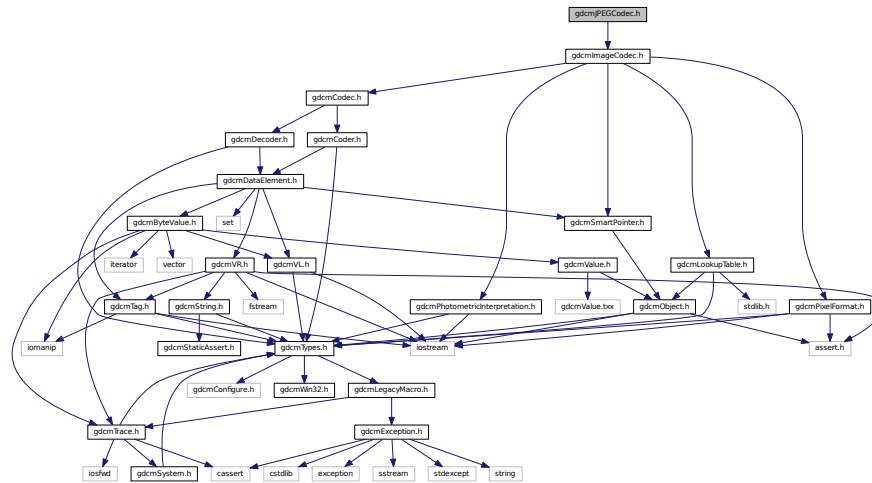

This graph shows which files directly or indirectly include this file:

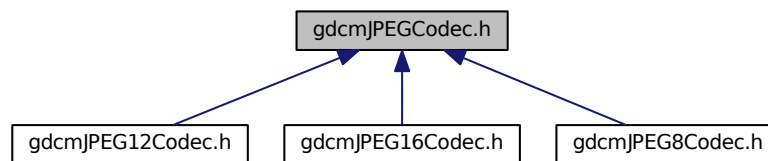

## Classes

- class [gdcm::JPEGCodec](#)

*JPEG codec Class to do JPEG (8bits, 12bits, 16bits lossy & lossless). It redispach in between the different codec implementation: [JPEG8Codec](#), [JPEG12Codec](#) & [JPEG16Codec](#) It also support inconsistency in between DICOM header and JPEG compressed stream [ImageCodec](#) implementation for the JPEG case.*

## Namespaces

- [gdcm](#)

## 28.137 gdcmJPEGLSCodec.h File Reference


```
#include "gdcmImageCodec.h"
```

[illegible]

- class `gdcm::JPEGLSCodec`  
*JPEG-LS.*

- **gdcm**

```
#include "gdcmFile.h"
#include "gdcmDataElement.h"
Include dependency graph for gdcmJSON.h:
```

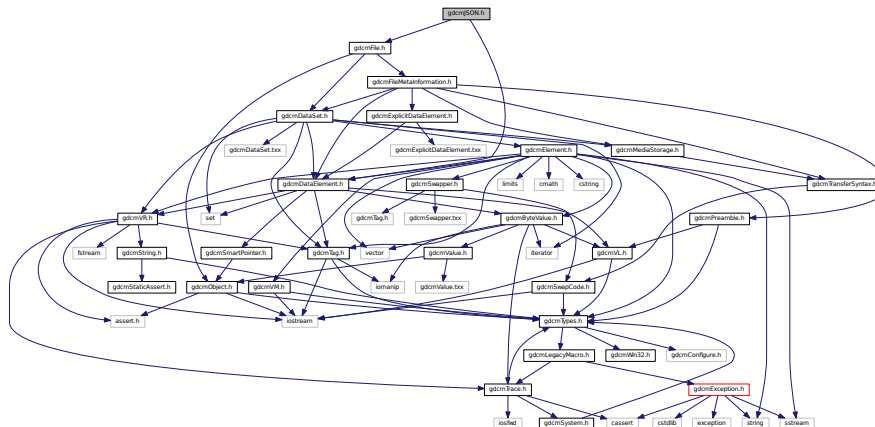


- class `gdc::JSON`

- **gdcm**

```
#include "gdcmImageCodec.h"
```

Include dependency graph for gdcMkKADUCodec.h:

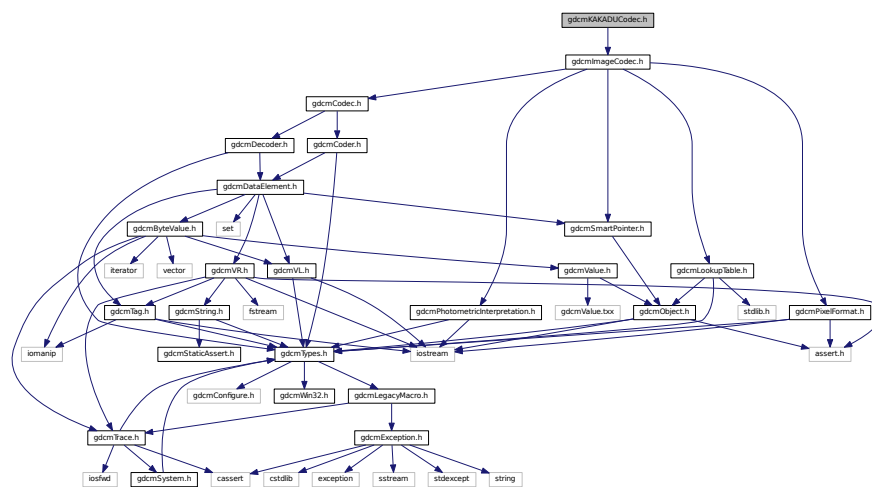


- class `gdcm::KAKADUCodec`  
*KAKADUCodec.*

- **gdcm**

```
#include "gdcmException.h"
```

```
#include "gdcmTrace.h"
```

Include dependency graph for gdcmLegacyMacro.h:

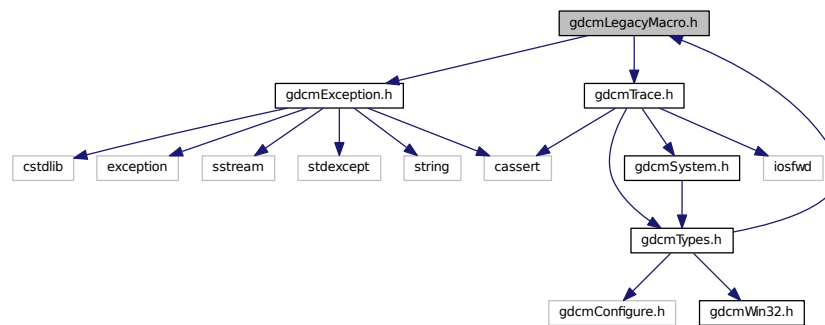

This graph shows which files directly or indirectly include this file:



## Macros

- `#define GDCM_LEGACY(method) method;`
- `#define GDCM_LEGACY_BODY(method, version) gdcmWarningMacro(#method " was deprecated for " version " and will be removed in a future version.")`
- `#define GDCM_LEGACY_REPLACED_BODY(method, version, replace) gdcmWarningMacro(#method " was deprecated for " version " and will be removed in a future version. Use " #replace " instead.")`

### 28.140.1 Macro Definition Documentation

28.140.1.1 `#define GDCM_LEGACY( method ) method;`


28.140.1.2 `#define GDCM_LEGACY_BODY( method, version ) gdcmWarningMacro(#method " was deprecated for " version " and will be removed in a future version.")`

28.140.1.3 `#define GDCM_LEGACY_REPLACED_BODY( method, version, replace ) gdcmWarningMacro(#method " was deprecated for " version " and will be removed in a future version. Use " #replace " instead.")`

## 28.141 gdcmLO.h File Reference

```
#include "gdcmString.h"
```

Include dependency graph for `gdcmLO.h`:

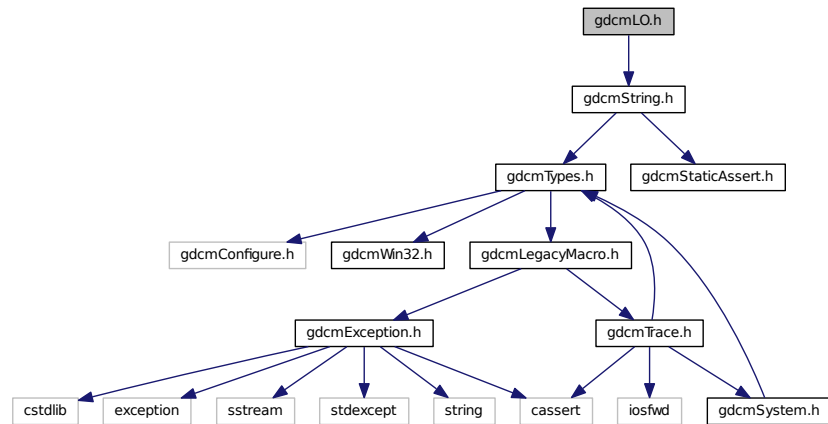

## Classes

- class `gdcm::LO`

*LO.*

## Namespaces

- `gdcm`


## 28.142 gdcmLookupTable.h File Reference

```


#include "gdcmTypes.h"
#include "gdcmObject.h"
#include <stdlib.h>

```

Include dependency graph for gdcmLookupTable.h:

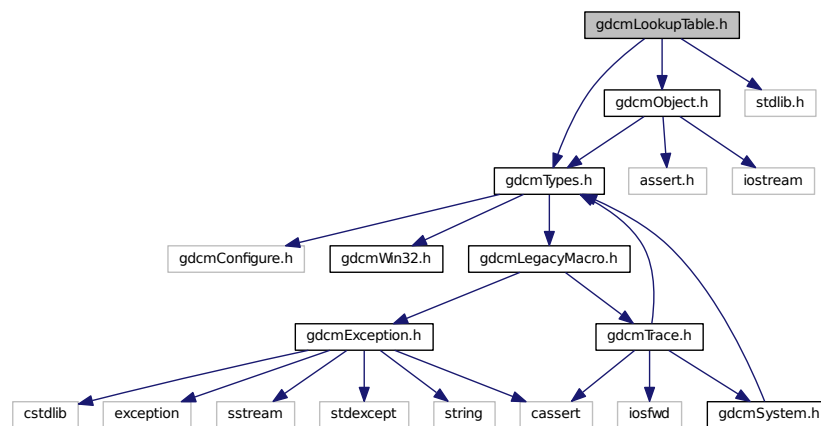

This graph shows which files directly or indirectly include this file:

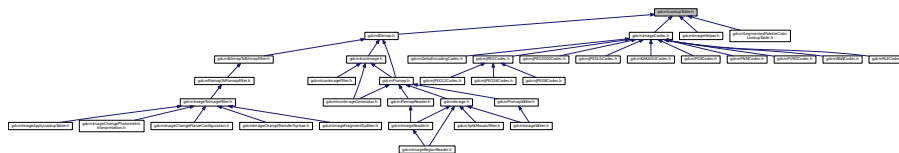

## Classes

- class [gdcm::LookupTable](#)  
*LookupTable* class.

## Namespaces

- [gdcm](#)


## 28.143 gdcmMacro.h File Reference

```


#include "gdcmTypes.h"
#include "gdcmTag.h"
#include "gdcmMacroEntry.h"
#include <map>
#include <vector>

```

Include dependency graph for gdcMacro.h:

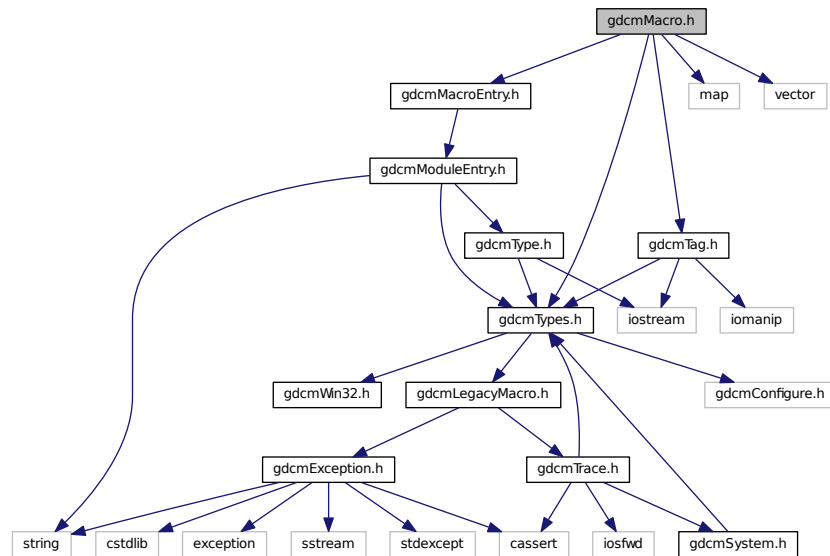

This graph shows which files directly or indirectly include this file:

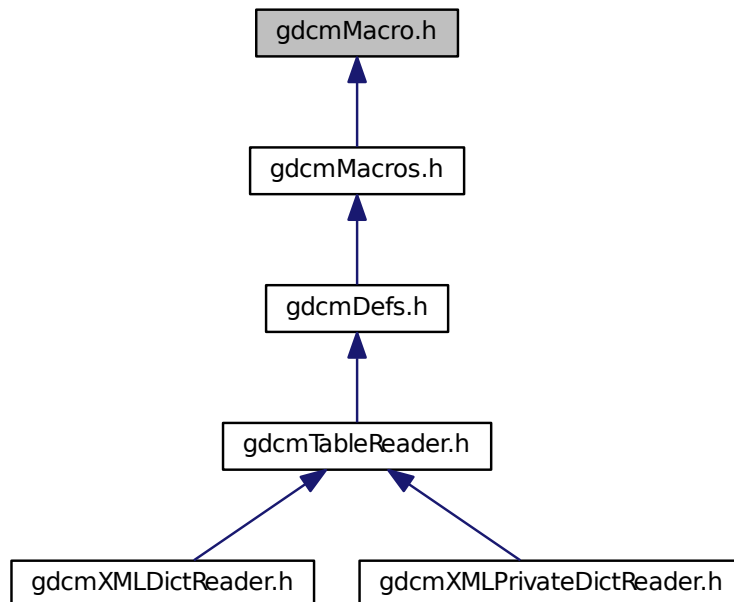



## Classes

- class [gdcm::Macro](#)

*Class for representing a [Macro](#).*

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const Macro &_val)`


## 28.144 gdcmMacroEntry.h File Reference

```
#include "gdcmModuleEntry.h"
```

Include dependency graph for gdcmMacroEntry.h:

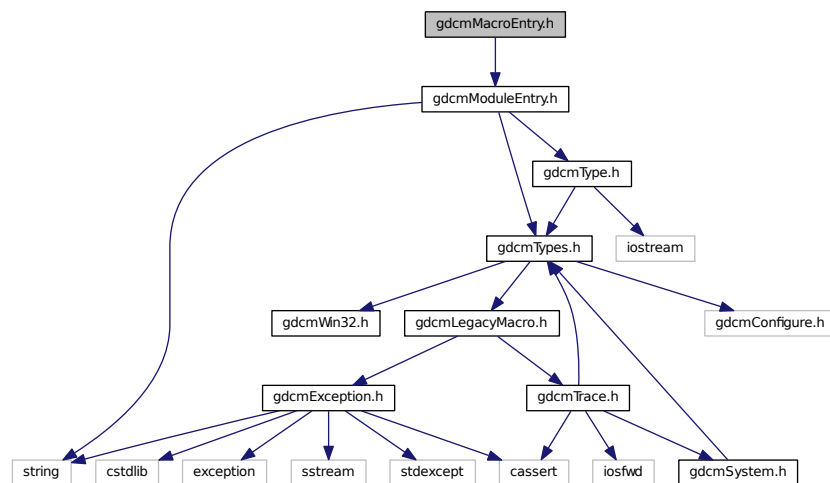

This graph shows which files directly or indirectly include this file:

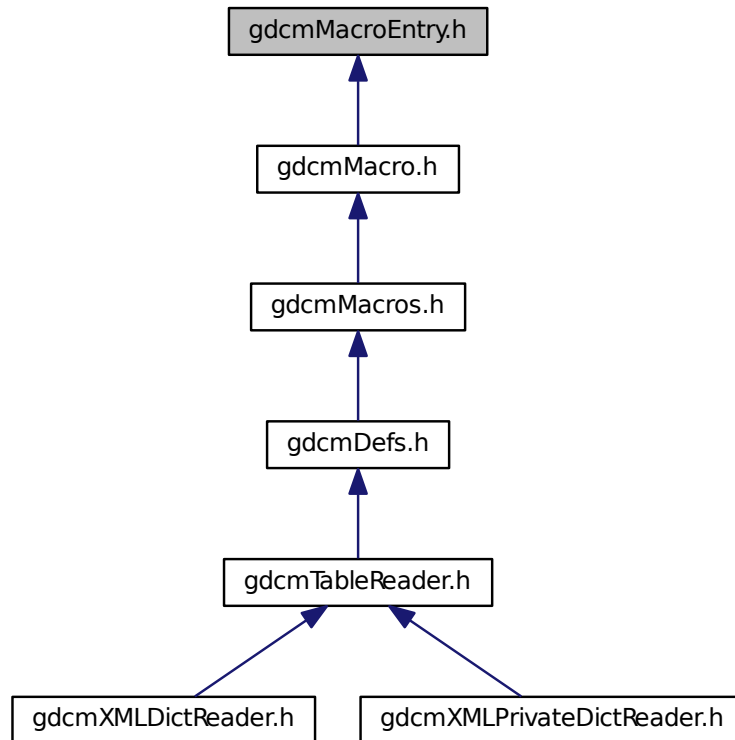

## Macros

- `#define` [GDCMMACROENTRY\\_H](#)


### 28.144.1 Macro Definition Documentation

#### 28.144.1.1 `#define` GDCMMACROENTRY\_H


## 28.145 gdcMacros.h File Reference

```
#include "gdcTypes.h"
#include "gdcMacro.h"
#include <map>
```

Include dependency graph for gdcmMacros.h:

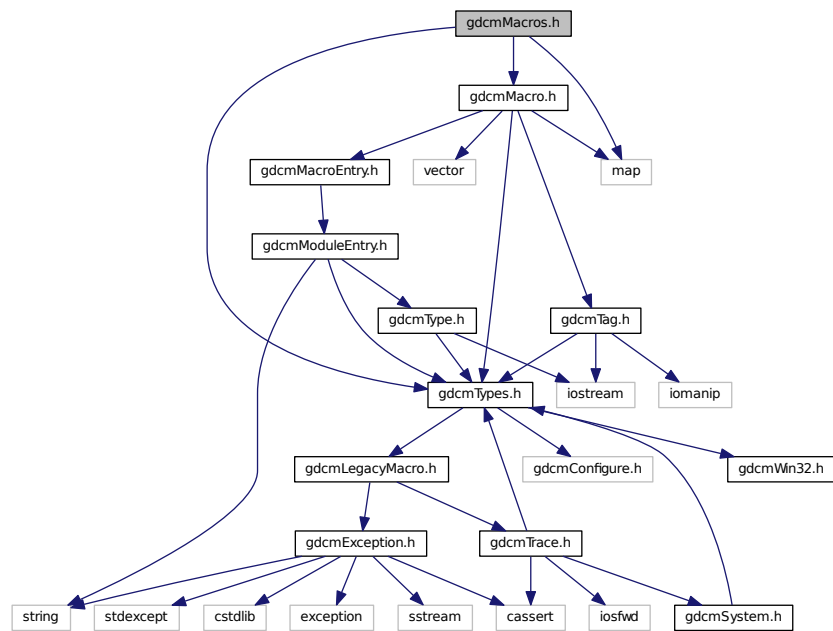

This graph shows which files directly or indirectly include this file:

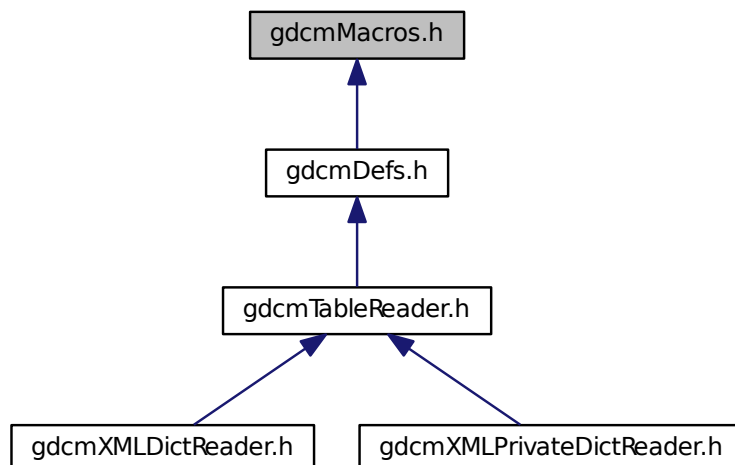

## Classes

- class [gdcm::Macros](#)

*Class for representing a [Modules](#).*

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const Macros &_val)`


## 28.146 gdcmMaximumLengthSub.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmMaximumLengthSub.h:

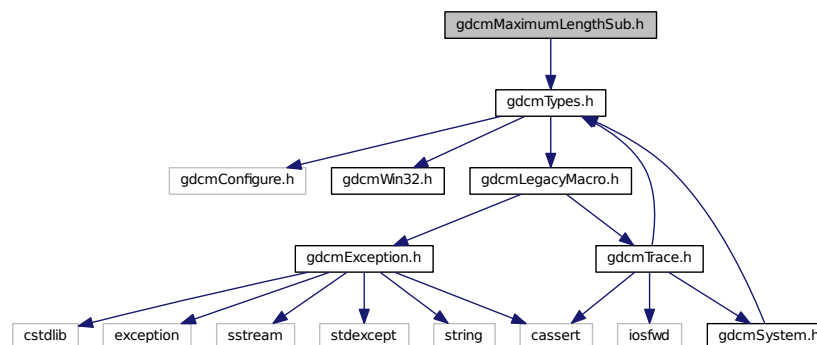

This graph shows which files directly or indirectly include this file:

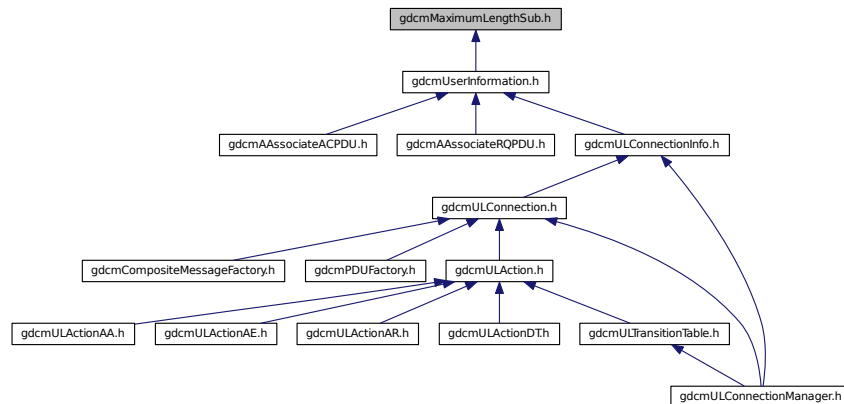

## Classes

- class [gdcm::network::MaximumLengthSub](#)  
[MaximumLengthSub](#) Annex D Table D.1-1 MAXIMUM LENGTH SUB-ITEM FIELDS (A-ASSOCIATE-RQ)


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.147 gdcmMD5.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmMD5.h:

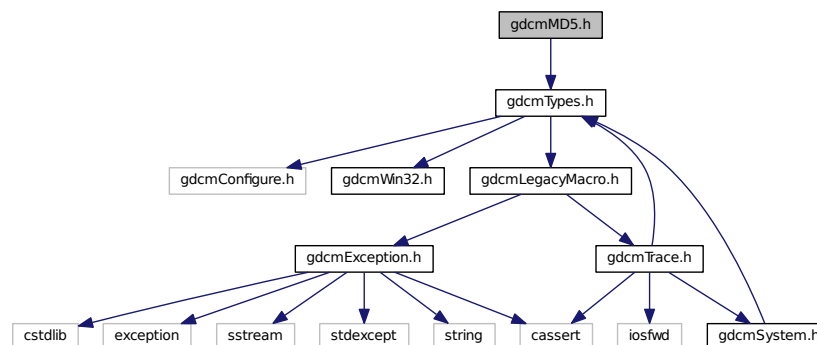

## Classes

- class [gdcm::MD5](#)  
Class for [MD5](#).


## Namespaces

- [gdcm](#)


## 28.148 gdcmMediaStorage.h File Reference

```
#include "gdcmTransferSyntax.h"
```

Include dependency graph for gdcmMediaStorage.h:

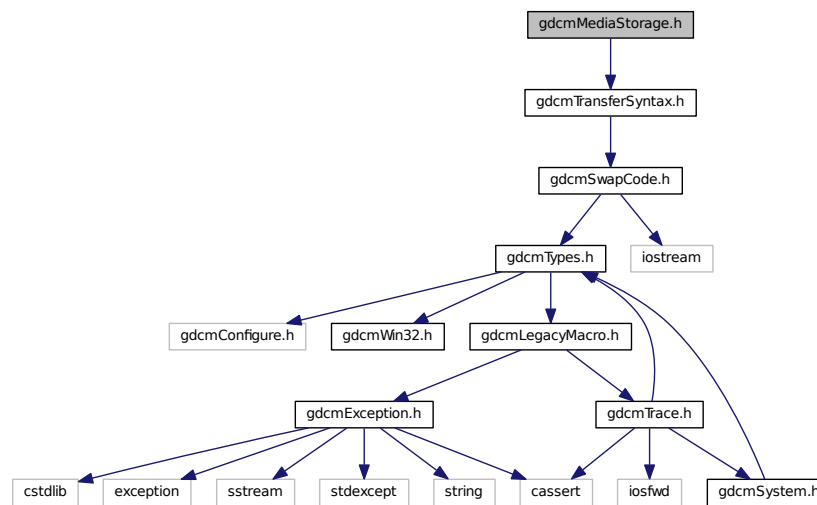

This graph shows which files directly or indirectly include this file:

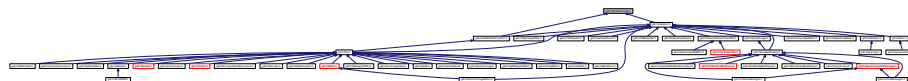

## Classes

- class [gdcm::MediaStorage](#)  
[MediaStorage](#).

## Namespaces


- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &_os, const MediaStorage &ms)`

## 28.149 gdcmMeshPrimitive.h File Reference

```
#include <gdcmObject.h>
#include <gdcmDataElement.h>
Include dependency graph for gdcmMeshPrimitive.h:
```

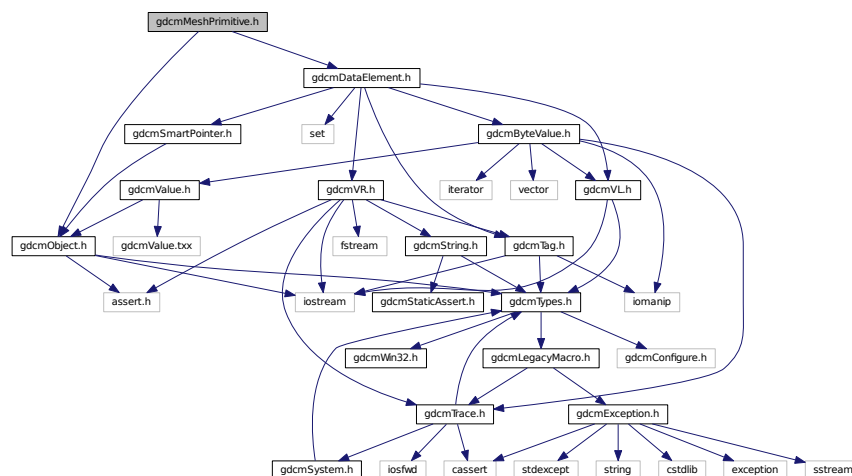

This graph shows which files directly or indirectly include this file:

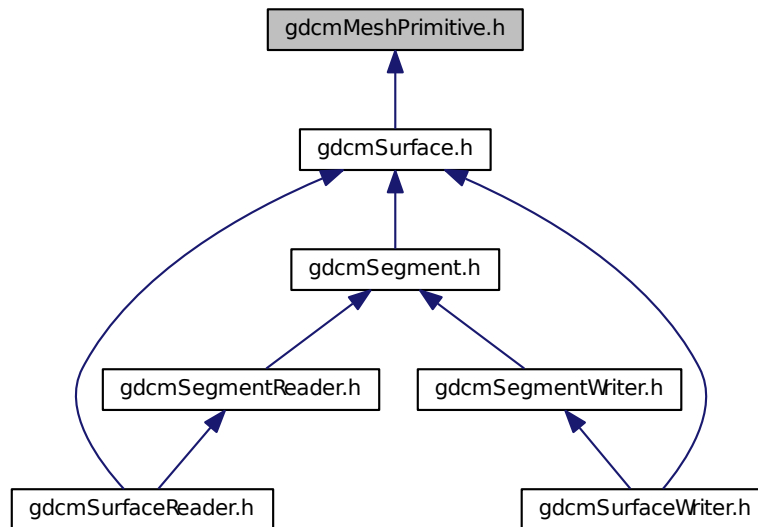

## Classes

- class [gdcm::MeshPrimitive](#)

*This class defines surface mesh primitives. It is designed from surface mesh primitives macro.*

## Namespaces

- [gdcm](#)

## 28.150 gdcmModule.h File Reference


```

#include "gdcmTypes.h"
#include "gdcmTag.h"
#include "gdcmModuleEntry.h"
#include <map>
#include <vector>


```



Include dependency graph for gdcmModule.h:

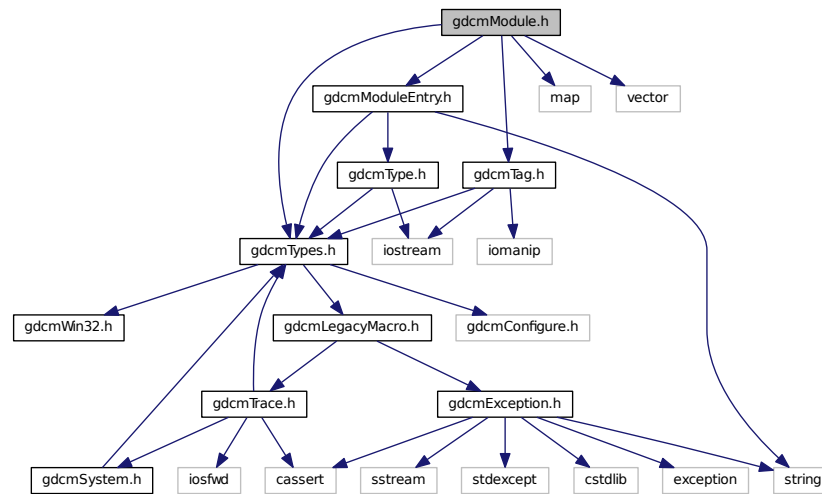


This graph shows which files directly or indirectly include this file:

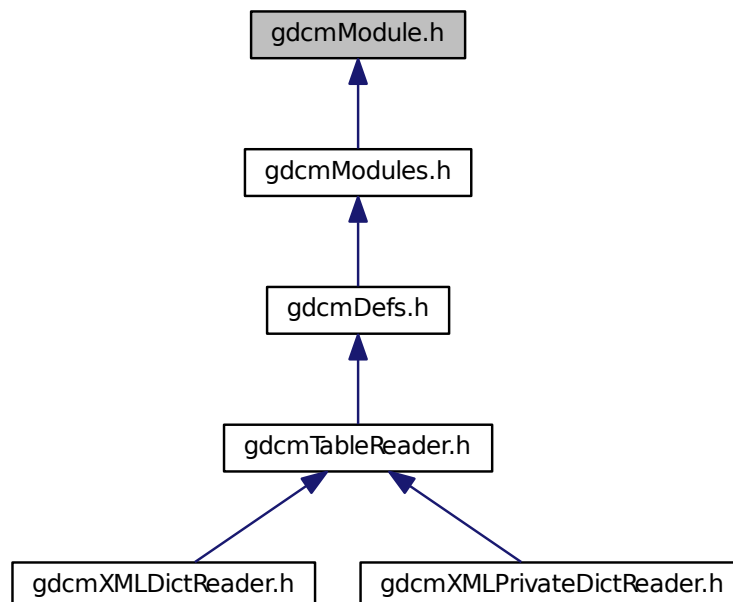



This graph shows which files directly or indirectly include this file:

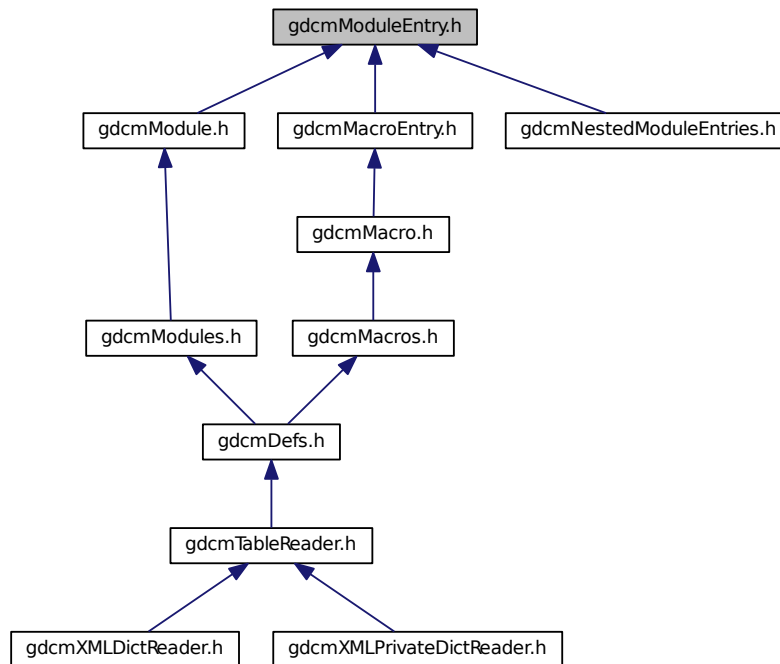

## Classes

- class [gdcm::ModuleEntry](#)  
Class for representing a *ModuleEntry*.

## Namespaces

- [gdcm](#)

## Typedefs

- typedef ModuleEntry [gdcm::MacroEntry](#)

## Functions

- std::ostream & [gdcm::operator<<](#) (std::ostream &\_os, const ModuleEntry &\_val)


## 28.152 gdcmModules.h File Reference

```
#include "gdcmTypes.h"
```


```
#include "gdcmModule.h"
```

```
#include <map>
```

Include dependency graph for gdcmModules.h:

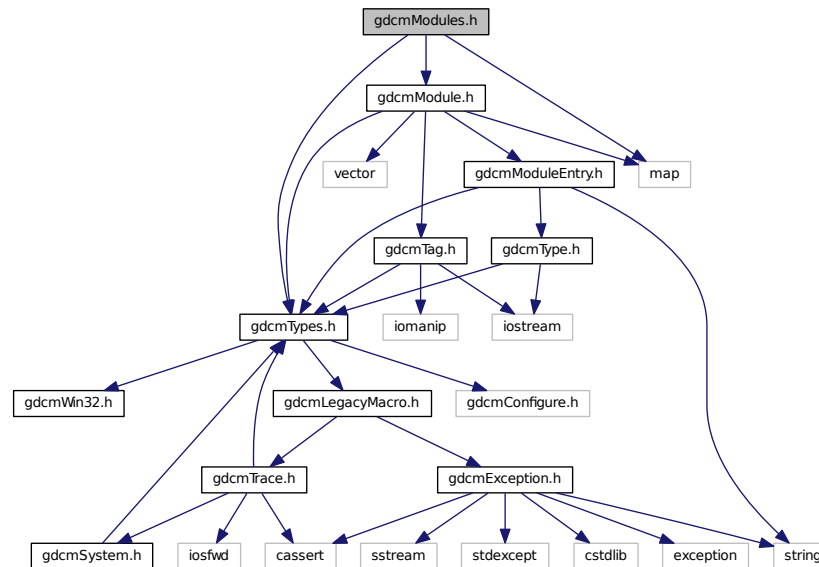


This graph shows which files directly or indirectly include this file:

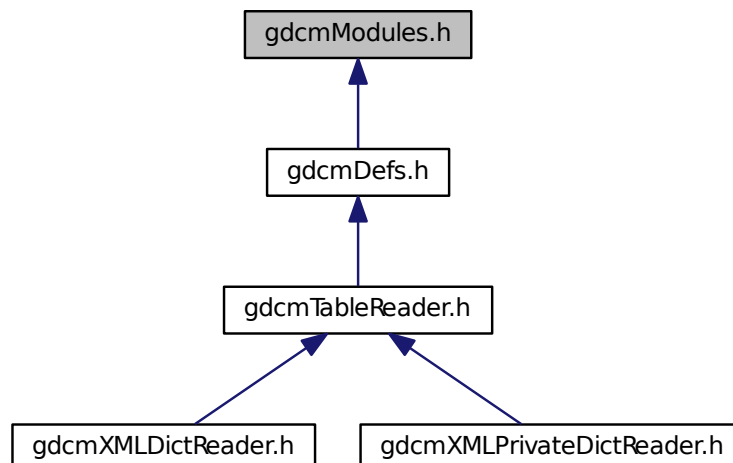





Include dependency graph for gdcmNestedModuleEntries.h:

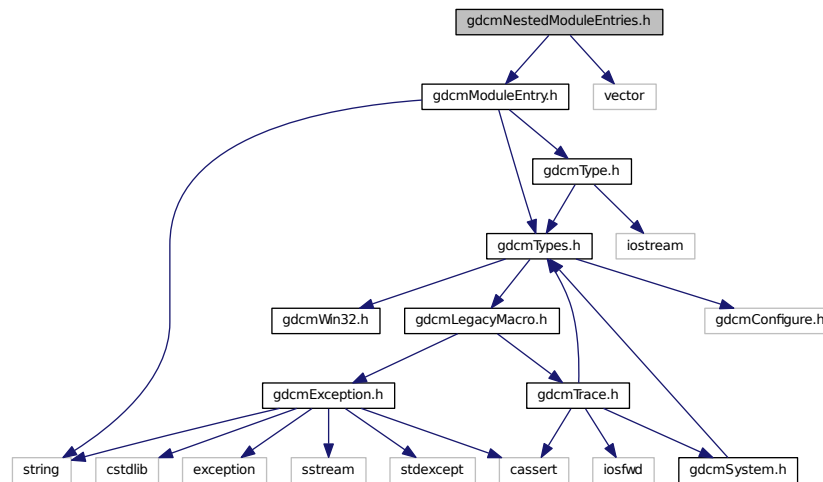

## Classes

- class [gdcm::NestedModuleEntries](#)

*Class for representing a [NestedModuleEntries](#).*

## Namespaces

- [gdcm](#)

## Typedefs


- typedef NestedModuleEntries [gdcm::NestedMacroEntries](#)

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const NestedModuleEntries &_val)`

## 28.156 gdcNetworkEvents.h File Reference

This graph shows which files directly or indirectly include this file:

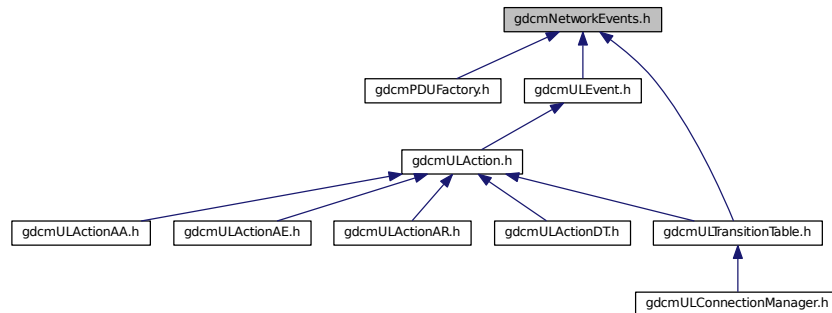

### Namespaces

- [gdc](#)
- [gdc::network](#)

### Enumerations

- `enum gdc::network::EEventID {`  
`gdc::network::eAASSOCIATERequestLocalUser = 0,`  
`gdc::network::eTransportConnConfirmLocal,`  
`gdc::network::eASSOCIATE_ACPDUreceived,`  
`gdc::network::eASSOCIATE_RJPDUreceived,`  
`gdc::network::eTransportConnIndicLocal,`  
`gdc::network::eAASSOCIATE_RQPDUreceived,`  
`gdc::network::eAASSOCIATEresponseAccept,`  
`gdc::network::eAASSOCIATEresponseReject,`  
`gdc::network::ePDATArequest,`  
`gdc::network::ePDATATFPDU,`  
`gdc::network::eARELEASERequest,`  
`gdc::network::eARELEASE_RQPDUReceivedOpen,`  
`gdc::network::eARELEASE_RPPDUReceived,`  
`gdc::network::eARELEASEResponse,`  
`gdc::network::eAABORTRequest,`  
`gdc::network::eAABORTPDUReceivedOpen,`  
`gdc::network::eTransportConnectionClosed,`  
`gdc::network::eARTIMTimerExpired,`  
`gdc::network::eUnrecognizedPDUReceived,`  
`gdc::network::eEventDoesNotExist }`


### Variables

- `const int gdc::network::cMaxEventID = eEventDoesNotExist`



## 28.157 gdcNetworkStateID.h File Reference

This graph shows which files directly or indirectly include this file:

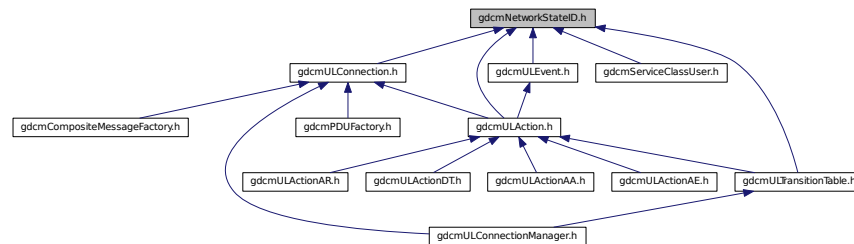

### Namespaces

- [gdc](#)
- [gdc::network](#)

### Enumerations

- enum [gdc::network::EStateID](#) {  
[gdc::network::eStaDoesNotExist](#) = 0,  
[gdc::network::eSta1Idle](#) = 1,  
[gdc::network::eSta2Open](#) = 2,  
[gdc::network::eSta3WaitLocalAssoc](#) = 4,  
[gdc::network::eSta4LocalAssocDone](#) = 8,  
[gdc::network::eSta5WaitRemoteAssoc](#) = 16,  
[gdc::network::eSta6TransferReady](#) = 32,  
[gdc::network::eSta7WaitRelease](#) = 64,  
[gdc::network::eSta8WaitLocalRelease](#) = 128,  
[gdc::network::eSta9ReleaseCollisionRqLocal](#) = 256,  
[gdc::network::eSta10ReleaseCollisionAc](#) = 512,  
[gdc::network::eSta11ReleaseCollisionRq](#) = 1024,  
[gdc::network::eSta12ReleaseCollisionAcLocal](#) = 2048,  
[gdc::network::eSta13AwaitingClose](#) = 4096 }

### Functions


- int [gdc::network::GetStateIndex](#) (EStateID inState)

### Variables


- const int [gdc::network::cMaxStateID](#) = 13

## 28.158 gdcmObject.h File Reference

```
#include "gdcmTypes.h"
#include <assert.h>
#include <iostream>
Include dependency graph for gdcmObject.h:
```

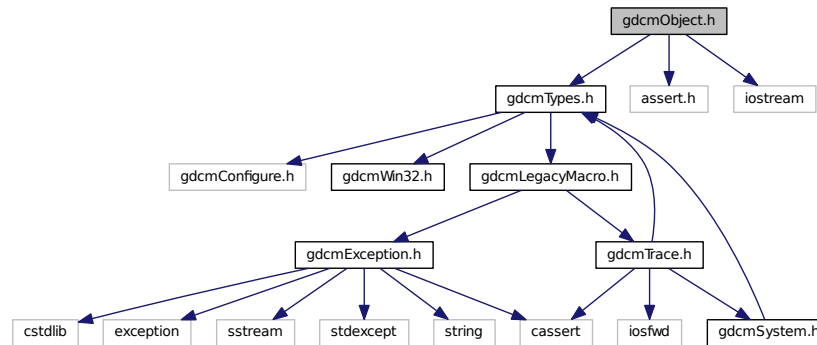

This graph shows which files directly or indirectly include this file:

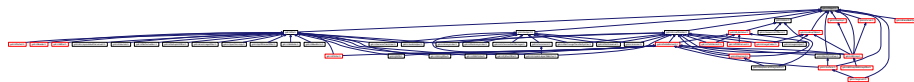

### Classes

- class [gdcm::Object](#)  
*Object.*
- class [gdcm::SmartPointer< ObjectType >](#)  
*Class for Smart Pointer.*

### Namespaces

- [gdcm](#)


### Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Object &obj)`

## 28.159 gdcmOpenSSLCryptoFactory.h File Reference

```
#include "gdcmCryptoFactory.h"
#include "gdcmOpenSSLCryptographicMessageSyntax.h"
```

Include dependency graph for gdcmOpenSSLCryptoFactory.h:

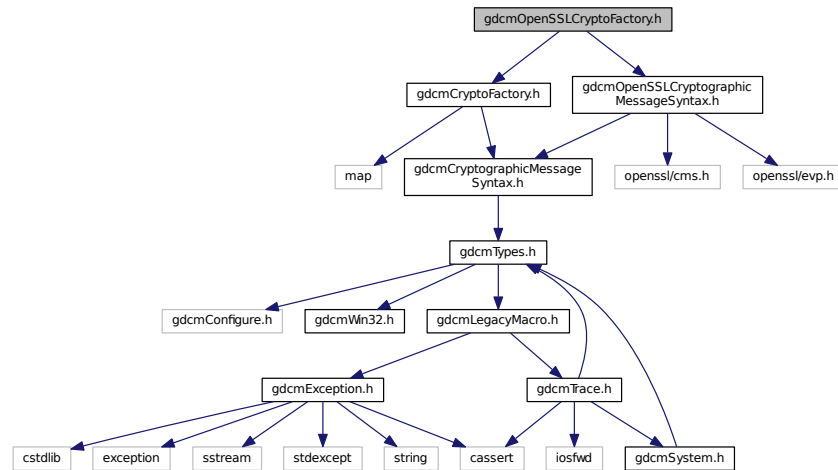

## Classes

- class [gdcm::OpenSSLCryptoFactory](#)

## Namespaces

- [gdcm](#)


## 28.160 gdcmOpenSSLCryptographicMessageSyntax.h File Reference

```


#include "gdcmCryptographicMessageSyntax.h"
#include <openssl/cms.h>
#include <openssl/evp.h>

```

Include dependency graph for `gdcmlOpenSSLCryptographicMessageSyntax.h`:

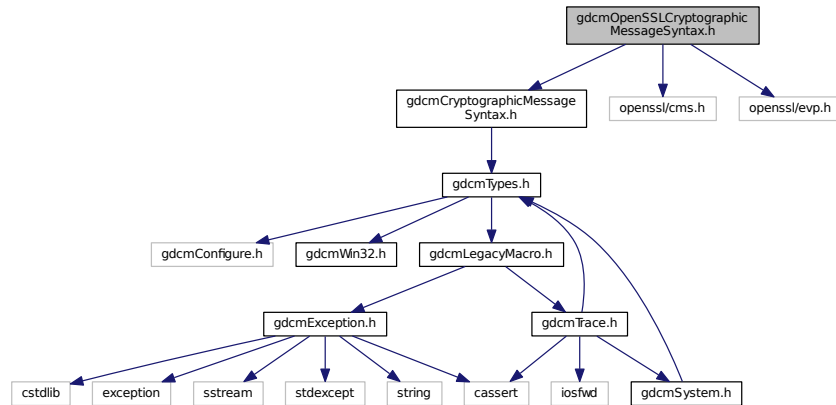

This graph shows which files directly or indirectly include this file:

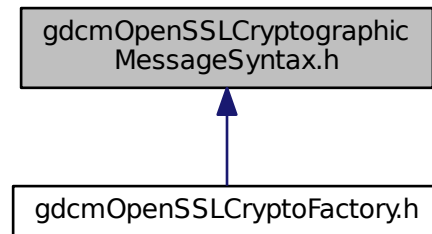

## Classes

- class [gdcml::OpenSSLCryptographicMessageSyntax](#)


## Namespaces

- [gdcml](#)

## 28.161 gdcmlOpenSSLP7CryptoFactory.h File Reference

```
#include "gdcmlCryptoFactory.h"
#include "gdcmlOpenSSLP7CryptographicMessageSyntax.h"
```

Include dependency graph for gdcmOpenSSLP7CryptoFactory.h:

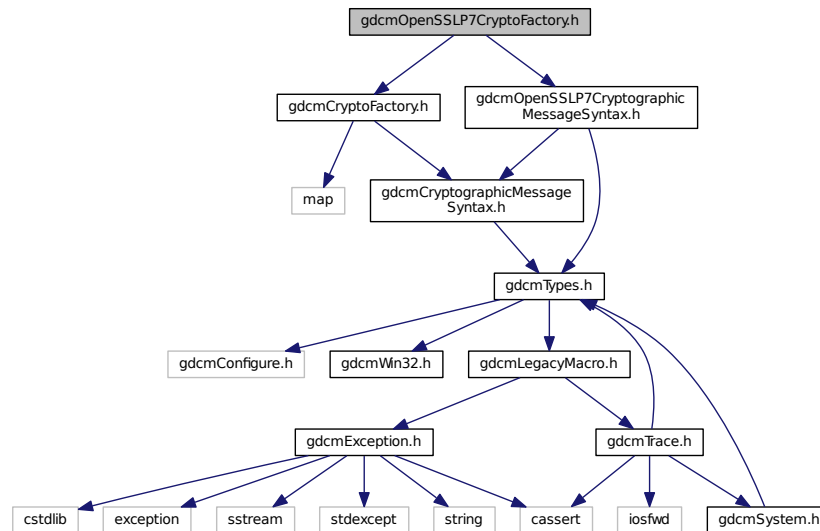

## Classes

- class [gdcm::OpenSSLP7CryptoFactory](#)


## Namespaces

- [gdcm](#)


## 28.162 gdcmOpenSSLP7CryptographicMessageSyntax.h File Reference

```
#include "gdcmCryptographicMessageSyntax.h"
#include "gdcmTypes.h"
```

Include dependency graph for `gdcOpenSSLP7CryptographicMessageSyntax.h`:

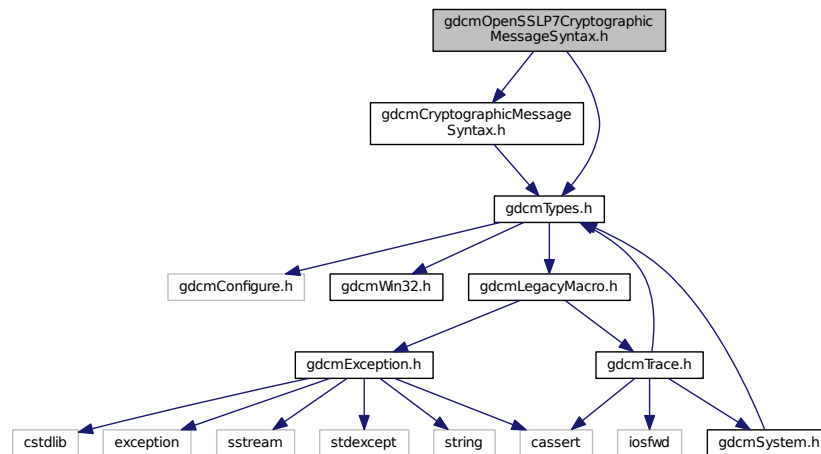

This graph shows which files directly or indirectly include this file:

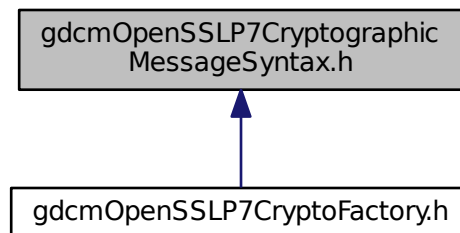

## Classes

- class [gdc::OpenSSLP7CryptographicMessageSyntax](#)

Class for [CryptographicMessageSyntax](#) encryption. This is just a simple wrapper around openssl PKCS7\_encrypt functionalities.


## Namespaces

- [gdc](#)

## 28.163 gdcmOrientation.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmOrientation.h:

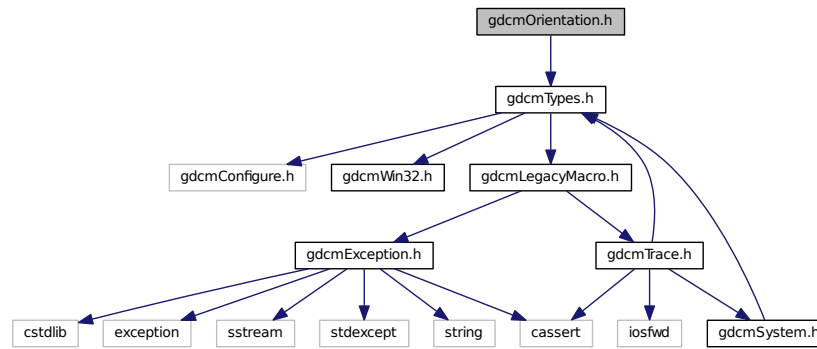

### Classes

- class [gdcm::Orientation](#)

*class to handle [Orientation](#)*

### Namespaces

- [gdcm](#)



### Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Orientation &o)`

## 28.164 gdcmOverlay.h File Reference

```
#include "gdcmTypes.h"
```

```
#include "gdcmObject.h"
```

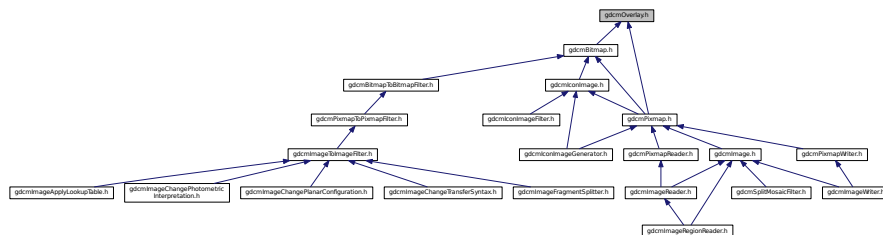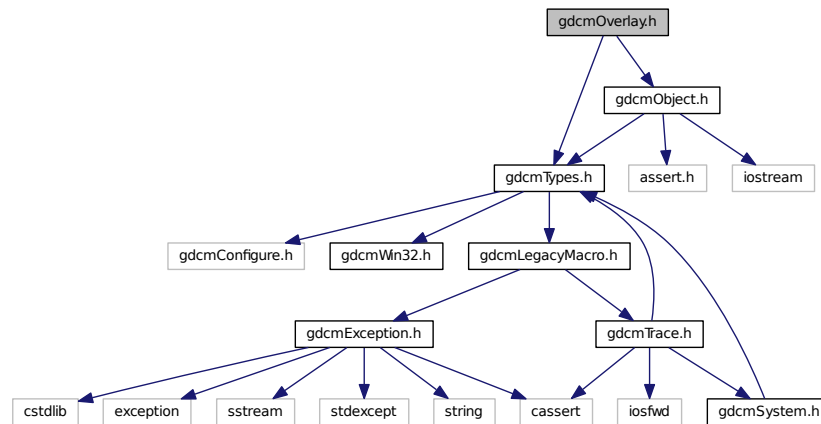

- class `gdcm::Overlay`  
*Overlay* class.


- **gdcm**

## 28.166 gdcmParseException.h File Reference


```
#include "gdcmException.h"
#include "gdcmDataElement.h"
```



Include dependency graph for gdcmParseException.h:

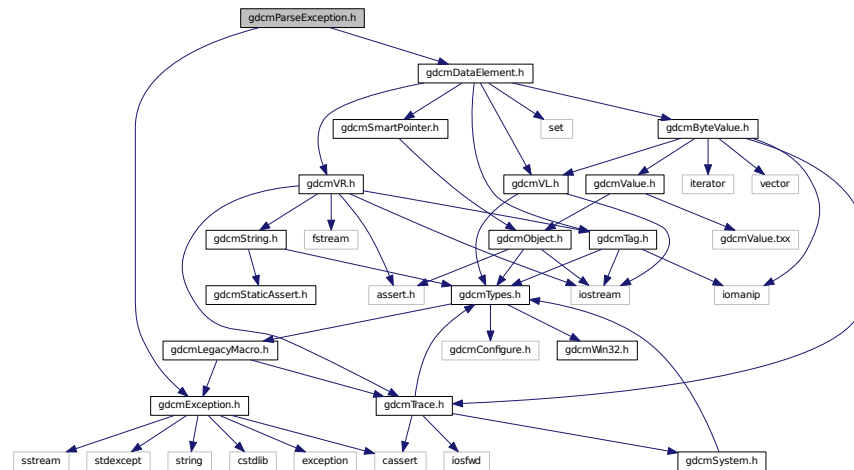

This graph shows which files directly or indirectly include this file:

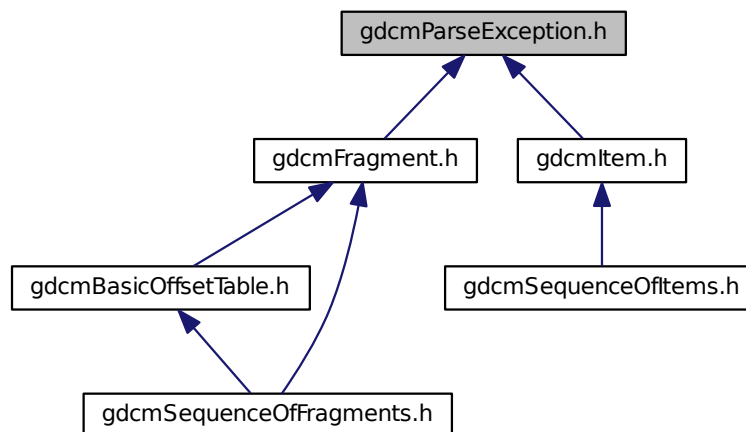

## Classes

- class [gdcm::ParseException](#)  
*ParseException* Standard exception handling object.


## Namespaces

- [gdcm](#)

## 28.167 gdcmParser.h File Reference

```
#include "gdcmTag.h"
#include "gdcmByteBuffer.h"
#include <fstream>
```

Include dependency graph for gdcmParser.h:

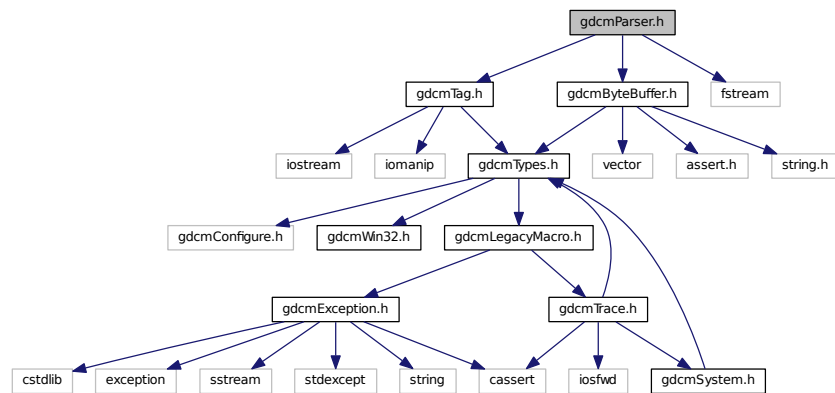

## Classes

- class [gdcm::Parser](#)

*Parser ala XML\_Parser from expat (SAX)*


## Namespaces

- [gdcm](#)

## 28.168 gdcmPatient.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmStudy.h"
```

Include dependency graph for gdcmPatient.h:

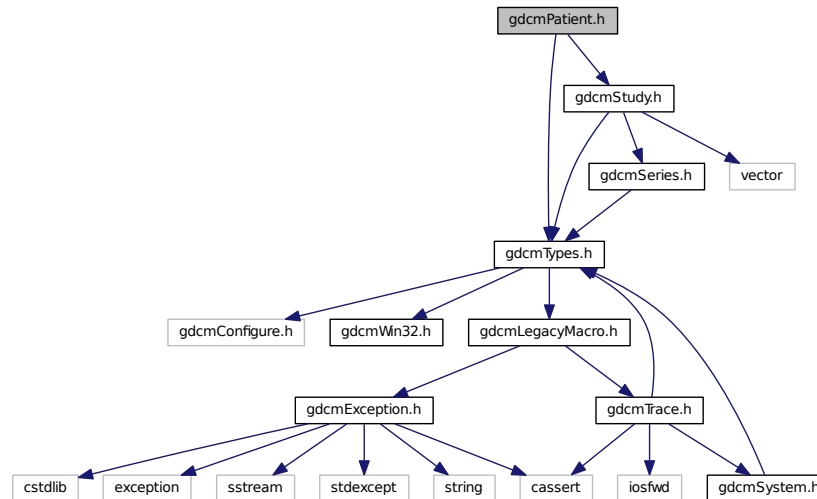

## Classes

- class [gdcm::Patient](#)

See PS 3.3 - 2007 DICOM MODEL OF THE REAL-WORLD, p 54.

## Namespaces

- [gdcm](#)


## 28.169 gdcmPDataTFPDU.h File Reference

```

#include "gdcmTypes.h"
#include "gdcmPresentationDataValue.h"
#include "gdcmBasePDU.h"
#include <limits>

```

Include dependency graph for `gdcmPidataTFPDU.h`:

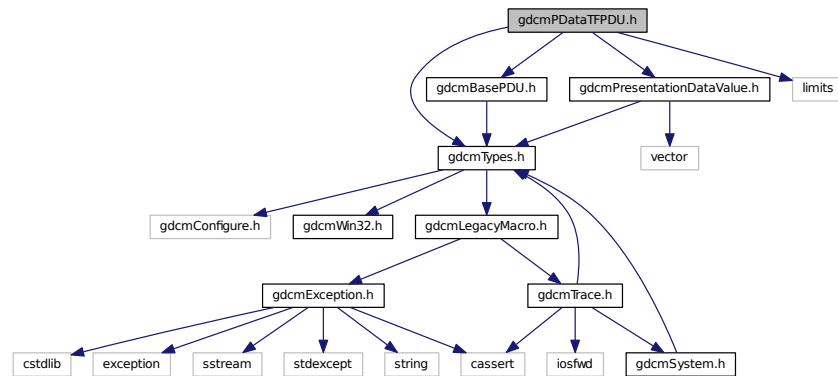

## Classes

- class `gdcm::network::PDataTFPDU`

*PDataTFPDU Table 9-22 P-DATA-TF PDU FIELDS.*

## Namespaces

- `gdcm`
- `gdcm::network`


## 28.170 gdcmPidBElement.h File Reference

```


#include "gdcmTag.h"
#include "gdcmVM.h"
#include "gdcmVR.h"
#include "gdcmByteValue.h"
#include "gdcmSmartPointer.h"

```

Include dependency graph for gdcmPDBElement.h:

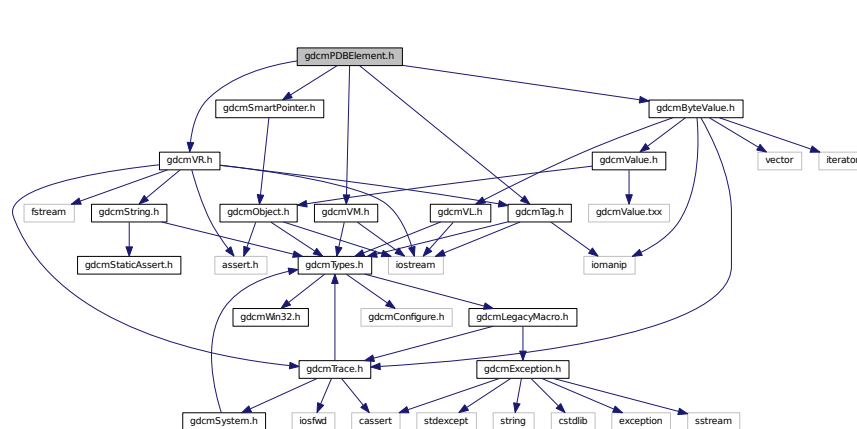

This graph shows which files directly or indirectly include this file:

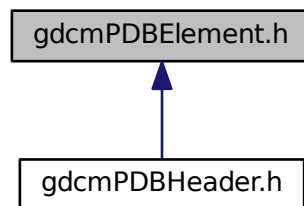

## Classes

- class [gdcm::PDBElement](#)  
*Class to represent a PDB [Element](#).*

## Namespaces


- [gdcm](#)

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const PDBElement &val)`



Include dependency graph for gdcmPDFCodec.h:

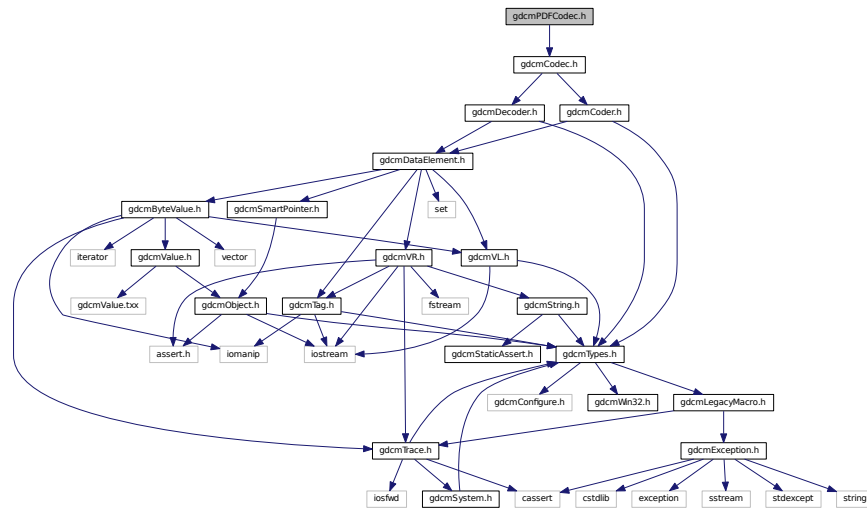

## Classes

- class [gdcm::PDFCodec](#)  
*PDFCodec* class.


## Namespaces

- [gdcm](#)

## 28.174 gdcmPDUFactory.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmNetworkEvents.h"
#include "gdcmULConnection.h"
#include "gdcmPresentationDataValue.h"
```

Include dependency graph for gdcmPDUFactory.h:

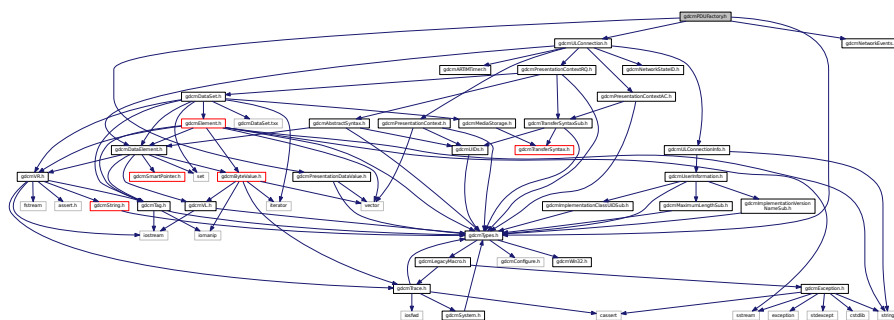

## Classes

- class [gdcm::network::PDUFactory](#)

*PDUFactory* basically, given an initial byte, construct the appropriate PDU. This way, the event loop doesn't have to know about all the different PDU types.


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.175 gdcmPersonName.h File Reference

```
#include "gdcmTypes.h"
#include <vector>
#include <algorithm>
#include <string.h>
```

Include dependency graph for gdcmPersonName.h:

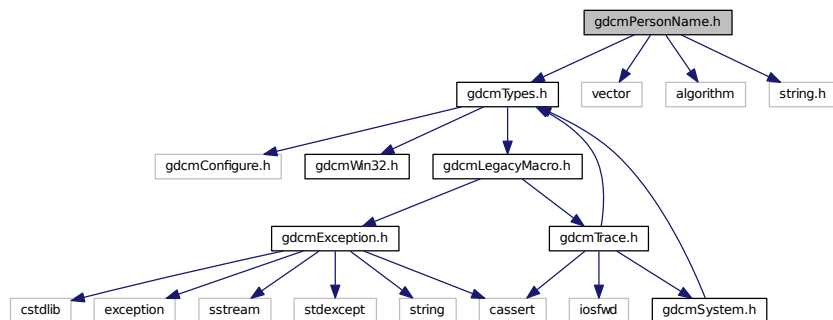

## Classes

- class [gdcm::PersonName](#)  
*PersonName* class.

## Namespaces

- [gdcm](#)

## 28.176 gdcmPGXCodec.h File Reference

```
#include "gdcmImageCodec.h"
```



- class `gdcm::PGXCodec`


- **gdcm**

```
#include "gdcmTypes.h"
#include <iostream>
```

```

graph TD
    gdcmPhotometricInterpretation.h --> gdcmTypes.h
    gdcmPhotometricInterpretation.h --> iostream
    gdcmTypes.h --> gdcmConfigure.h
    gdcmTypes.h --> gdcmWin32.h
    gdcmTypes.h --> gdcmLegacyMacro.h
    gdcmTypes.h --> gdcmException.h
    gdcmTypes.h --> gdcmTrace.h
    gdcmLegacyMacro.h --> gdcmException.h
    gdcmLegacyMacro.h --> gdcmTrace.h
    gdcmException.h --> cstdlb
    gdcmException.h --> exception
    gdcmException.h --> sstream
    gdcmException.h --> stdexcept
    gdcmException.h --> string
    gdcmException.h --> cassert
    gdcmException.h --> iosfwd
    gdcmException.h --> gdcmSystem.h
    gdcmTrace.h --> iosfwd
    gdcmTrace.h --> gdcmSystem.h
    gdcmSystem.h --> iosfwd
  
```

This graph shows which files directly or indirectly include this file:

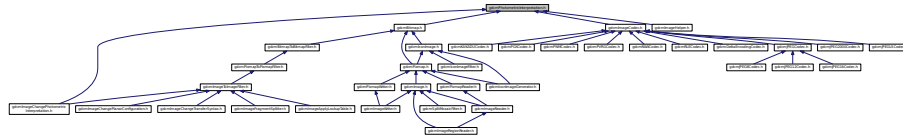

## Classes

- class [gdcm::PhotometricInterpretation](#)  
Class to represent an *PhotometricInterpretation*.

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const PhotometricInterpretation &val)`


## 28.178 gdcmPixelFormat.h File Reference

```
#include "gdcmTypes.h"
#include <iostream>
#include <assert.h>
```

Include dependency graph for gdcmPixelFormat.h:

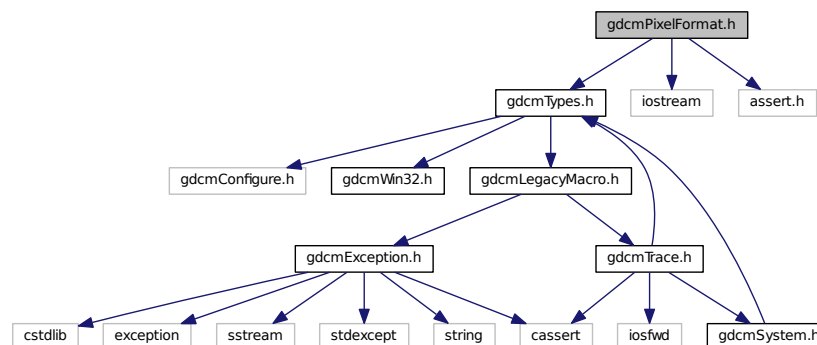

This graph shows which files directly or indirectly include this file:

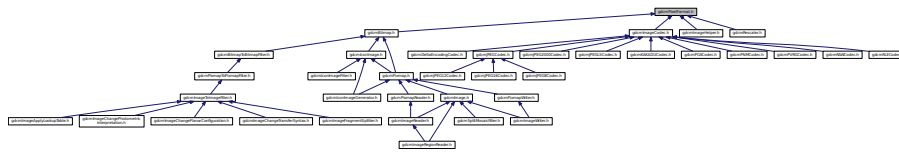

## Classes

- class [gdcm::PixelFormat](#)  
*PixelFormat.*

## Namespaces

- [gdcm](#)


## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const PixelFormat &pf)`

## 28.179 gdcmPixmap.h File Reference

```
#include "gdcmBitmap.h"
#include "gdcmCurve.h"
#include "gdcmIconImage.h"
#include "gdcmOverlay.h"
```

Include dependency graph for gdcmPixmap.h:

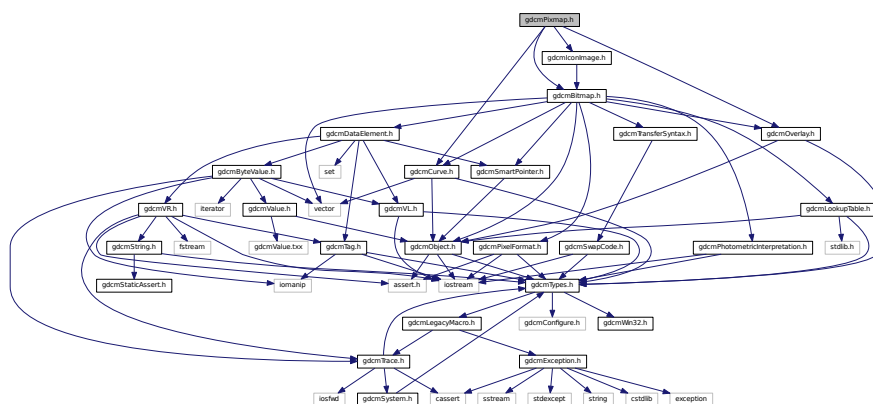

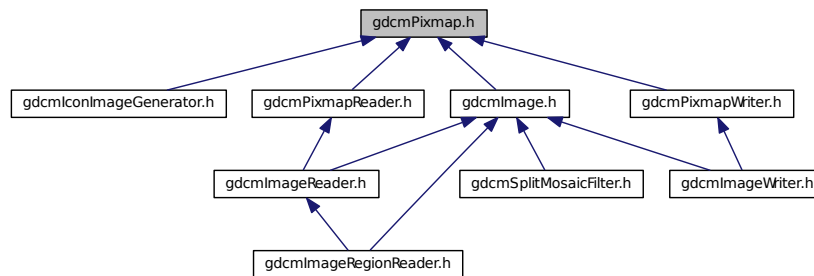


- class `gdcm::Pixmap`

*Bitmap class A bitmap based image. Used as parent for both IconImage and the main Pixel Data Image It does not contains any World Space information (IPP, IOP)*

- **gdcm**

```
#include "gdcmReader.h"
#include "gdcmPixmap.h"
```

Include dependency graph for gdcmPidxmapReader.h:

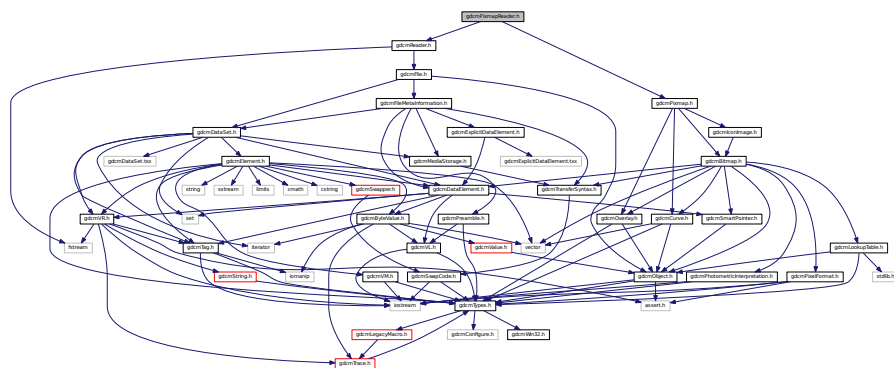



```
graph BT; gdcmlImageRegionReader.h --> gdcmlImageReader.h; gdcmlImageReader.h --> gdcmPixmapReader.h
```

- class `gdcm::PixmapReader`  
*PixmapReader*.

- **gdcm**


```
#include "gdcmBitmapToBitmapFilter.h"
Include dependency graph for gdcmPixmapToPixmapFilter.h:
```

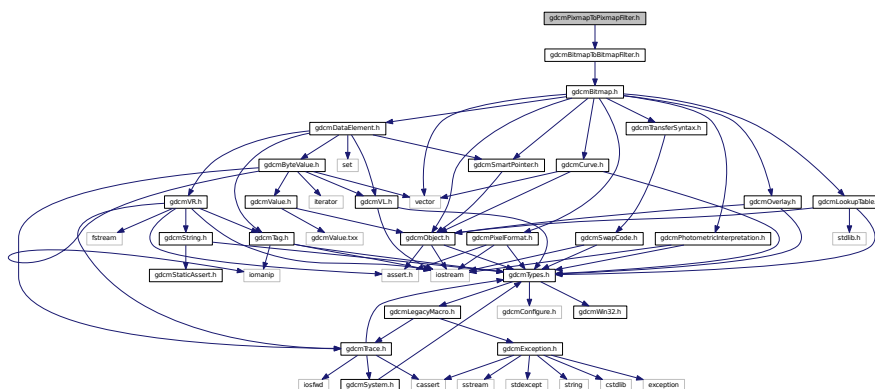

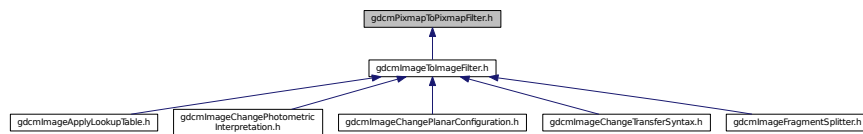

- PixmapToPixmapFilter* class Super class for all filter taking an image and producing an output image.

```
#include "gdcmPixmap.h"
```

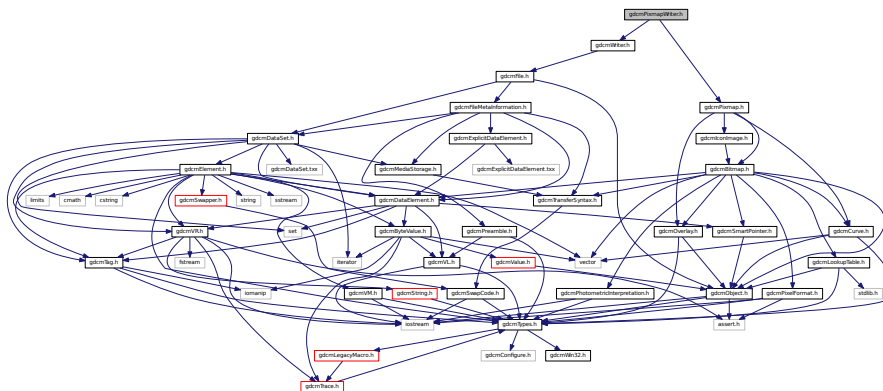


```
graph BT; gdcmImageWriter.h --> gdcmPixmapWriter.h
```

- class `gdcm::PixmapWriter`  
*PixmapWriter* This class will takes two inputs:

- `gdcm`

```
#include "gdcImageCodec.h"
Include dependency graph for gdcPNMCodec.h:
```

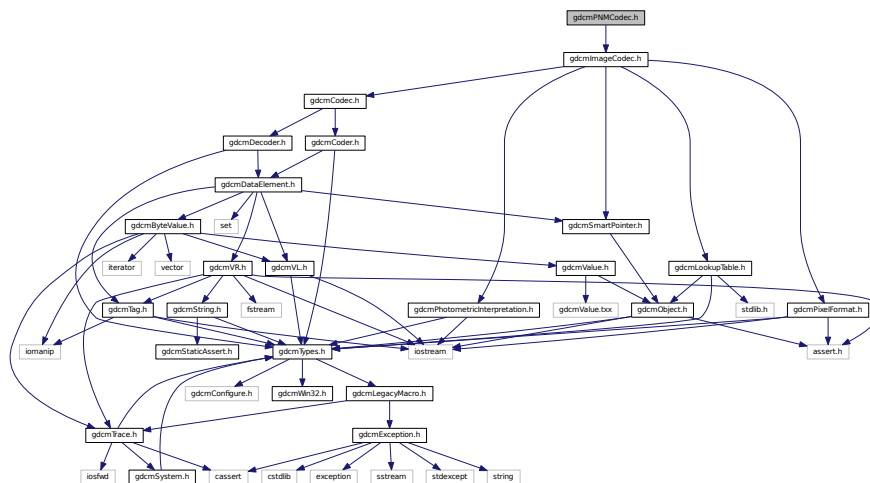

## Classes

- class [gdcm::PNMCodec](#)

*Class to do PNM PNM is the Portable anymap file format. The main web page can be found at: <http://netpbm.sourceforge.net/>.*

## Namespaces


- [gdcm](#)

## 28.184 gdcmPreamble.h File Reference


```
#include "gdcmTypes.h"
```

```
#include "gdcmVL.h"
```

Include dependency graph for gdcmPreamble.h:

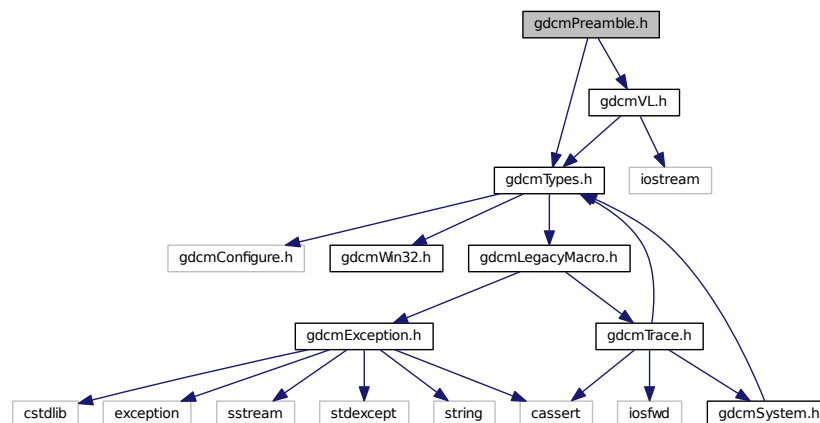

This graph shows which files directly or indirectly include this file:

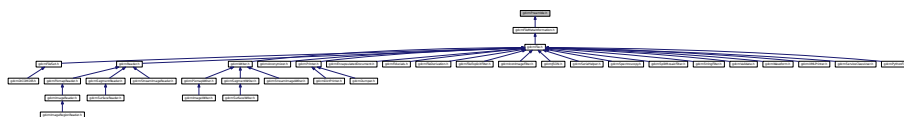

## Classes

- class [gdcm::Preamble](#)

*DICOM Preamble (Part 10)*



## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Preamble &val)`


## 28.185 gdcmPresentationContext.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmUIDs.h"
#include <vector>
```

Include dependency graph for gdcmPresentationContext.h:

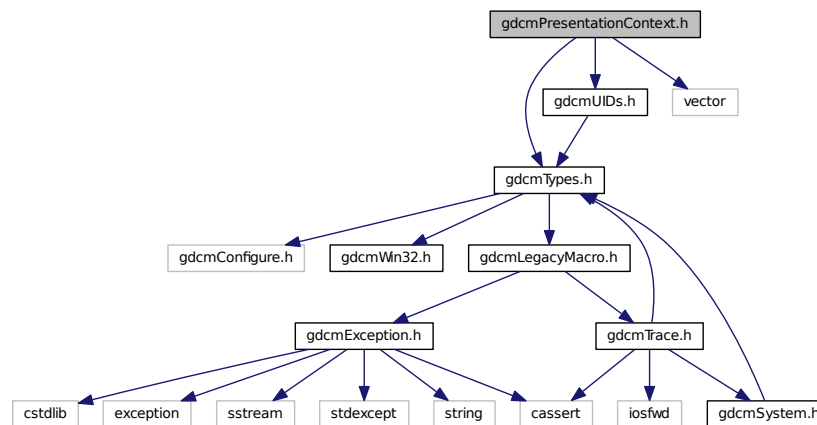

This graph shows which files directly or indirectly include this file:

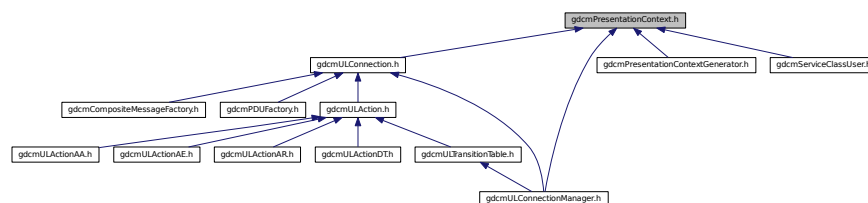

## Classes

- class [gdcm::PresentationContext](#)  
*PresentationContext.*

## Namespaces


- [gdcm](#)

## 28.186 gdcmPresentationContextAC.h File Reference


```
#include "gdcmTypes.h"
```

```
#include "gdcmTransferSyntaxSub.h"
```

Include dependency graph for gdcmPresentationContextAC.h:

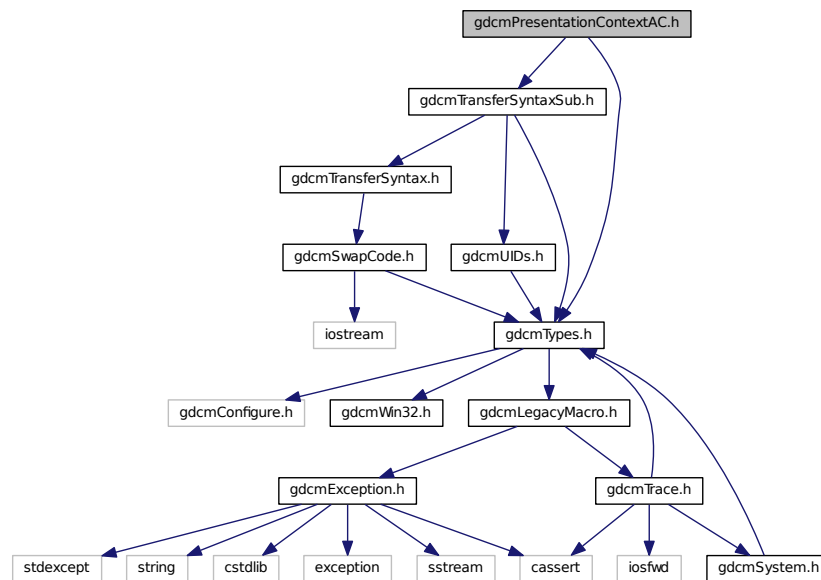

This graph shows which files directly or indirectly include this file:

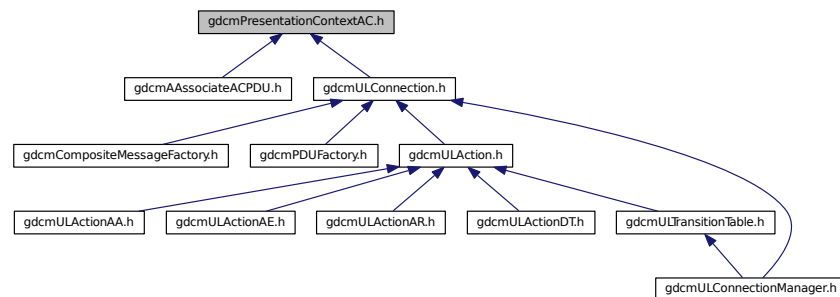

## Classes

- class [gdcm::network::PresentationContextAC](#)

[PresentationContextAC](#) [Table 9-18 PRESENTATION CONTEXT ITEM FIELDS.](#)

## Namespaces


- [gdcm](#)
- [gdcm::network](#)

## 28.187 gdcmPresentationContextGenerator.h File Reference

```
#include "gdcmDirectory.h"
```

```
#include "gdcmPresentationContext.h"
```

Include dependency graph for gdcmPresentationContextGenerator.h:

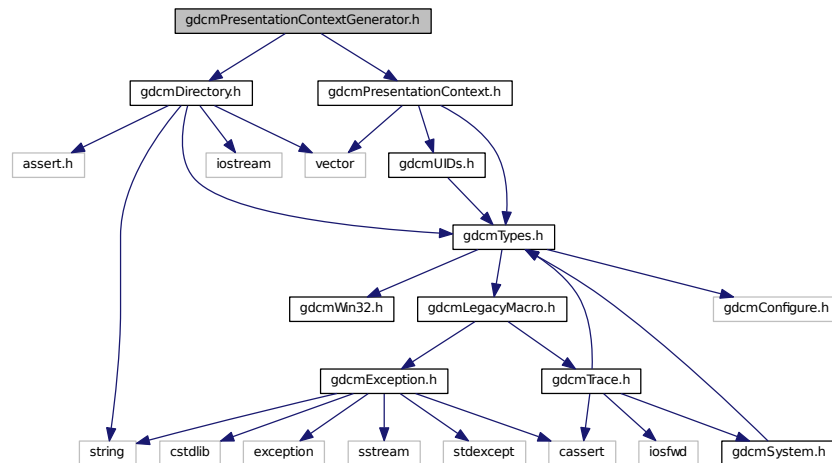

## Classes

- class [gdcm::PresentationContextGenerator](#)

[PresentationContextGenerator](#) This class is responsible for generating the proper [PresentationContext](#) that will be used in subsequent operation during a DICOM Query/Retrieve association. The step of the association is very sensible as special care need to be taken to explicitly define what instance are going to be send and how they are encoded.

## Namespaces

- [gdcm](#)




## 28.189 gdcmPresentationDataValue.h File Reference


```
#include "gdcmTypes.h"
```

```
#include <vector>
```

Include dependency graph for gdcmPresentationDataValue.h:

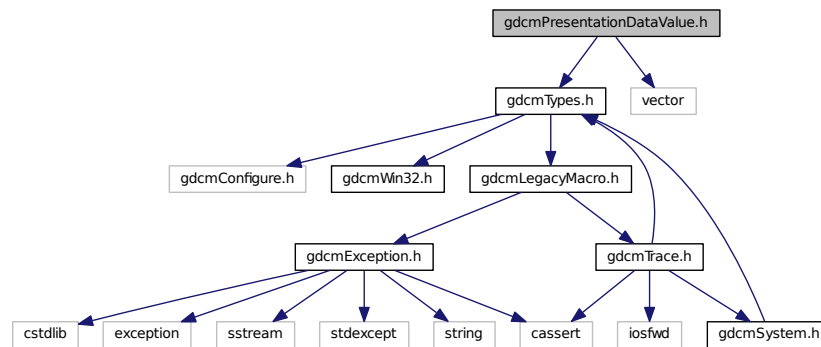

This graph shows which files directly or indirectly include this file:

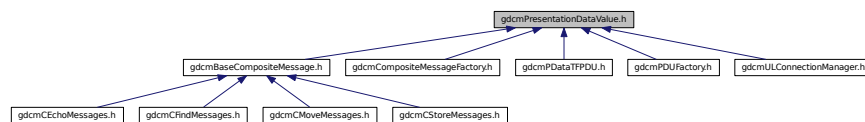

### Classes

- class [gdcm::network::PresentationDataValue](#)

*PresentationDataValue Table 9-23 PRESENTATION-DATA-VALUE ITEM FIELDS.*

### Namespaces


- [gdcm](#)
- [gdcm::network](#)

## 28.190 gdcmPrinter.h File Reference


```
#include "gdcmFile.h"
```

```
#include "gdcmDataElement.h"
```

Include dependency graph for `gdcmPrinter.h`:

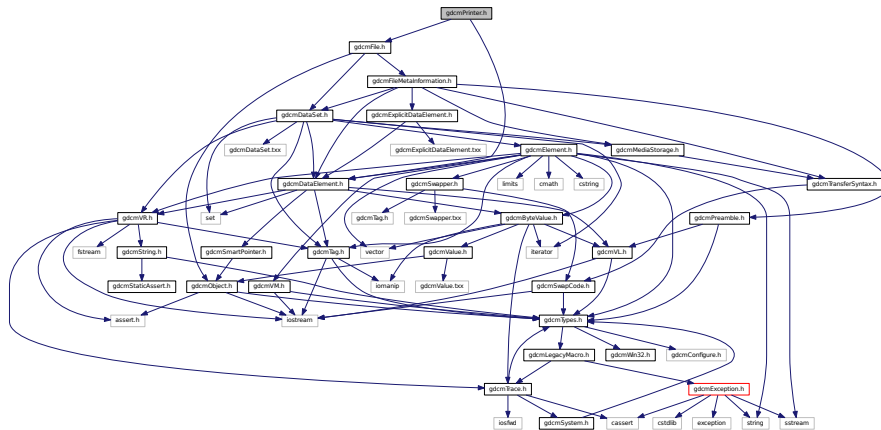

This graph shows which files directly or indirectly include this file:

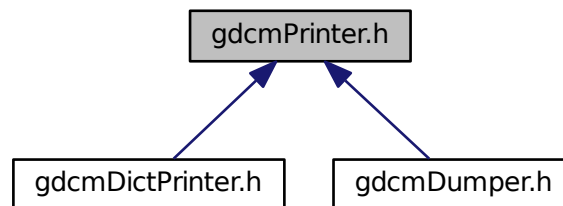

## Classes

- class `gdcm::Printer`  
*Printer* class.

## Namespaces

- **gdcm**

## 28.191 gdcmPrivateTag.h File Reference


```
#include "gdcmTag.h"
```

```


#include "gdcVR.h"
#include "gdcDataElement.h"
#include <iostream>
#include <iomanip>
#include <string>
#include <algorithm>
#include <string.h>
#include <ctype.h>

```

Include dependency graph for gdcPrivateTag.h:

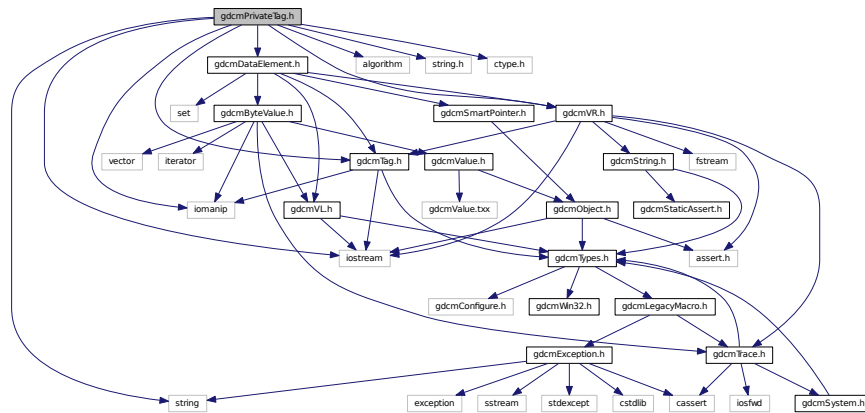

This graph shows which files directly or indirectly include this file:

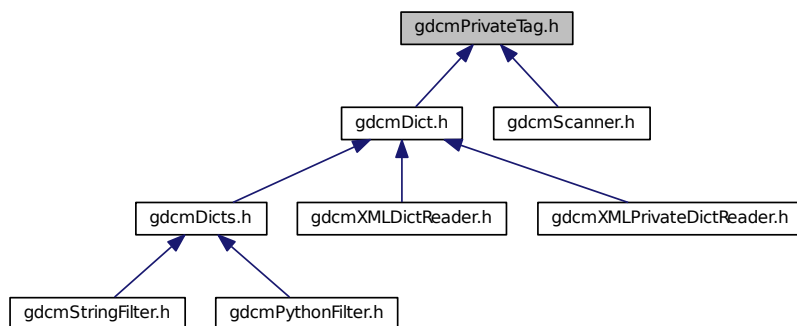

## Classes

- class [gdcm::PrivateTag](#)

Class to represent a Private DICOM Data [Element](#) ([Attribute](#)) [Tag](#) ([Group](#), [Element](#), [Owner](#))

## Namespaces

- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const PrivateTag &val)`

## 28.192 gdcmProgressEvent.h File Reference

```
#include "gdcmEvent.h"
```

```
#include "gdcmTag.h"
```

Include dependency graph for `gdcmProgressEvent.h`:

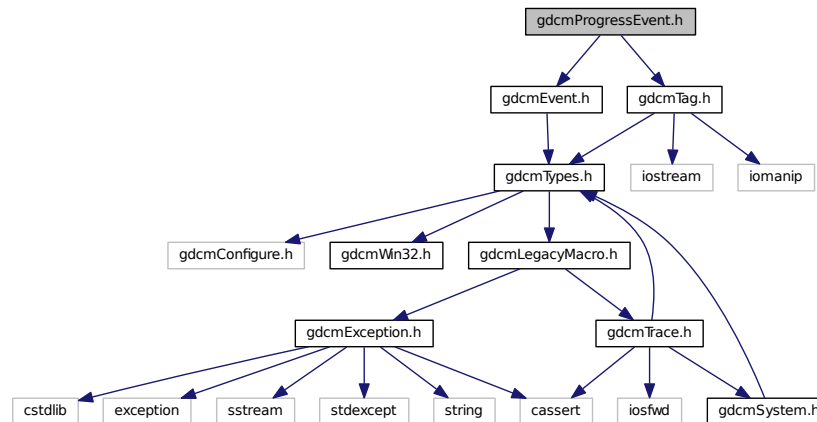

## Classes

- class [gdcm::ProgressEvent](#)  
*ProgressEvent* Special type of event triggered during.

## Namespaces

- [gdcm](#)

## 28.193 gdcmPVRGCodec.h File Reference

```
#include "gdcmImageCodec.h"
```



[illegible]

- class `gdcm::PVRGCodec`  
*PVRGCodec*.

- **gdcm**

```
#include <Python.h>
#include "gdcmDataElement.h"
#include "gdcmDicts.h"
#include "gdcmFile.h"
```

## Classes

- class [gdcm::PythonFilter](#)

*PythonFilter* [PythonFilter](#) is the class that make *gdcm2.x* looks more like *gdcm1* and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.

## Namespaces

- [gdcm](#)


## 28.195 gdcmQueryBase.h File Reference

```
#include "gdcmTag.h"
```


```
#include "gdcmDataElement.h"
```

```
#include <vector>
```

Include dependency graph for `gdcmQueryBase.h`:

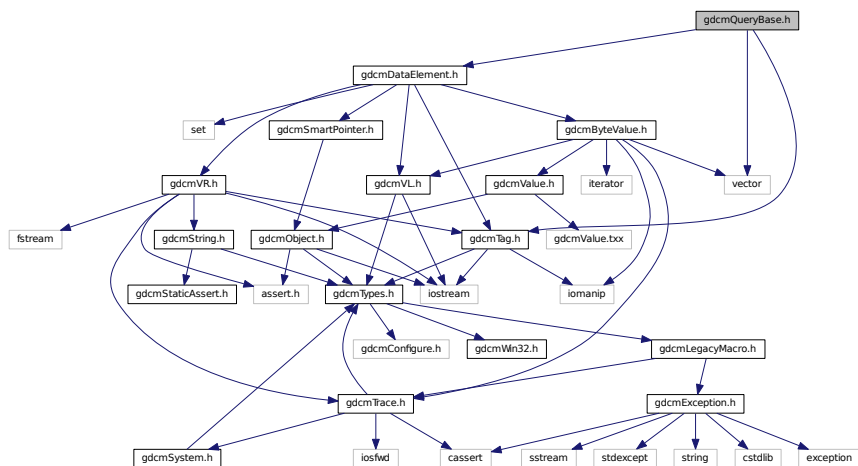

This graph shows which files directly or indirectly include this file:

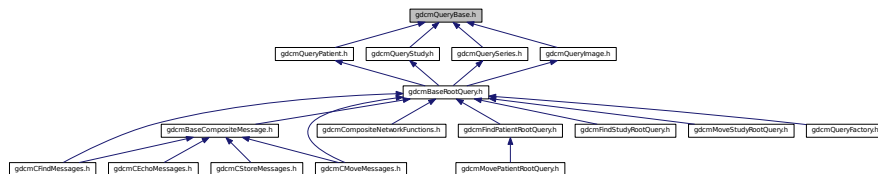

## Classes


- class [gdcm::QueryBase](#)

*QueryBase* contains: the base class for constructing a query dataset for a C-FIND and a C-MOVE.

- **gdcm**

- enum `gdcm::ERootType` {  
`gdcm::ePatientRootType`,  
`gdcm::eStudyRootType` }

```
#include "gdcmBaseRootQuery.h"
Include dependency graph for gdcmQueryFactory.h:
```

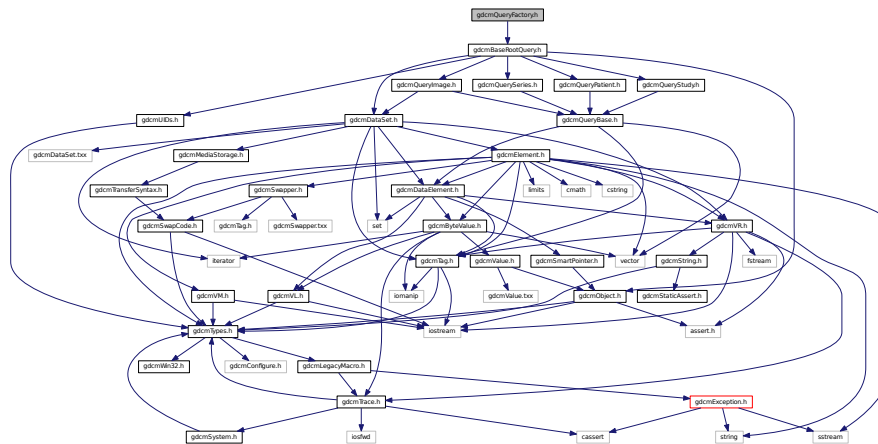

- class `gdcm::QueryFactory`  
*QueryFactory.h*

- **gdcm**

## Enumerations

```

• enum gdcmm::ECharSet {
    gdcmm::eLatin1 = 0,
    gdcmm::eLatin2,
    gdcmm::eLatin3,
    gdcmm::eLatin4,
    gdcmm::eCyrillic,
    gdcmm::eArabic,
    gdcmm::eGreek,
    gdcmm::eHebrew,
    gdcmm::eLatin5,
    gdcmm::eJapanese,
    gdcmm::eThai,
    gdcmm::eJapaneseKanjiMultibyte,
    gdcmm::eJapaneseSupplementaryKanjiMultibyte,
    gdcmm::eKoreanHangulHanjaMultibyte,
    gdcmm::eUTF8,
    gdcmm::eGB18030 }


```

## 28.197 gdcmmQueryImage.h File Reference


```
#include "gdcmmQueryBase.h"
```

```
#include "gdcmmDataSet.h"
```

Include dependency graph for gdcmmQueryImage.h:

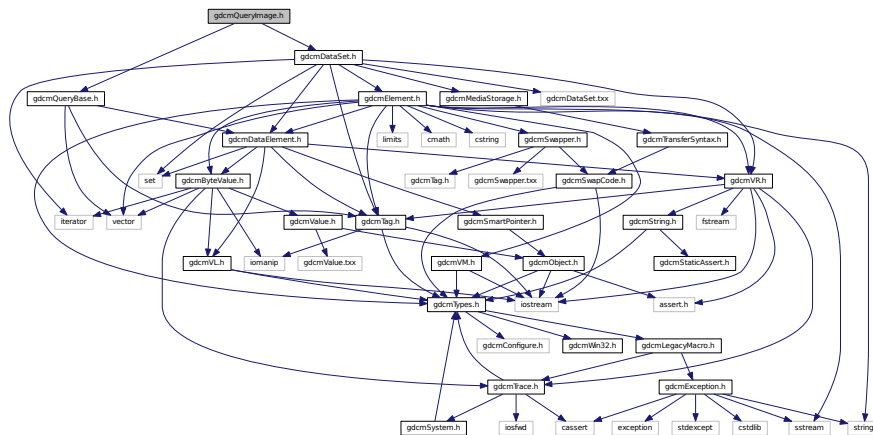

This graph shows which files directly or indirectly include this file:

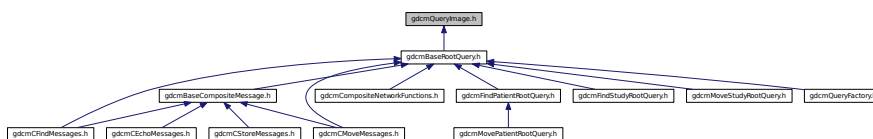

## Classes

- class `gdcm::QueryImage`

*QueryImage* contains: class to construct an image-based query for C-FIND and C-MOVE.


## Namespaces

- **gdcm**


## 28.198 gdcQueryPatient.h File Reference

```
#include "gdcmQueryBase.h"
```

Include dependency graph for gdcmQueryPatient.h:

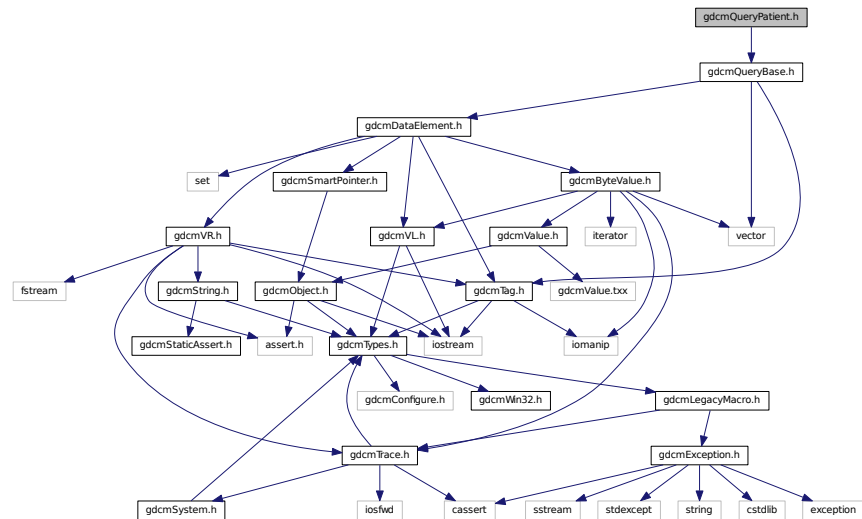

This graph shows which files directly or indirectly include this file:

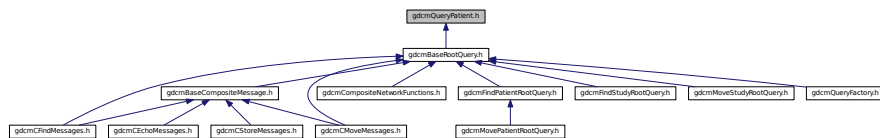

## Classes

- class `gdcm::QueryPatient`

*QueryPatient* contains: class to construct a patient-based query for c-find and c-move.


## Namespaces

- [gdcm](#)


## 28.199 gdcmQuerySeries.h File Reference

```
#include "gdcmQueryBase.h"
```

Include dependency graph for gdcmQuerySeries.h:

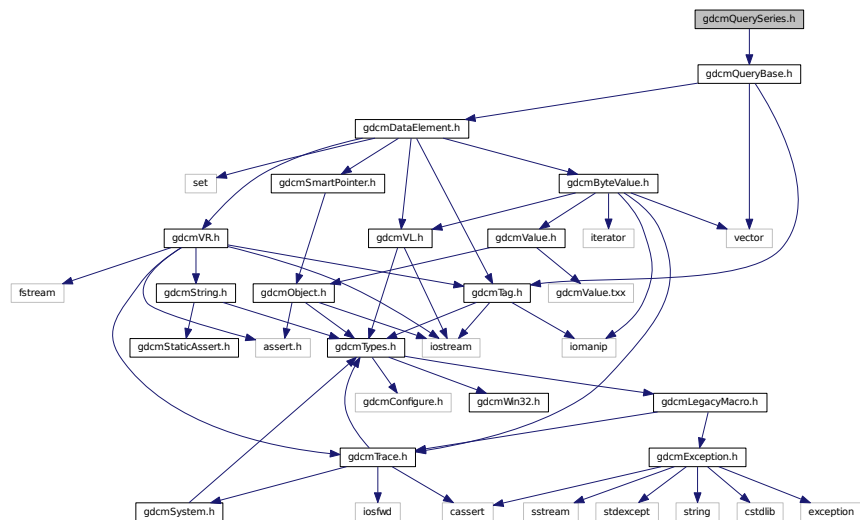

This graph shows which files directly or indirectly include this file:

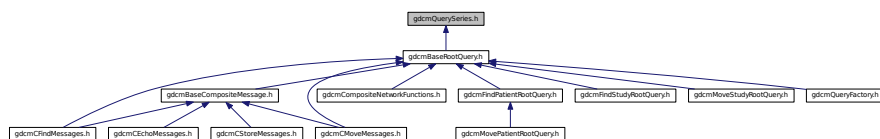

## Classes

- class [gdcm::QuerySeries](#)  
*QuerySeries* contains: class to construct a series-based query for c-find and c-move.


## Namespaces

- [gdcm](#)


## 28.200 gdcmQueryStudy.h File Reference

```
#include "gdcmQueryBase.h"
```

Include dependency graph for gdcmQueryStudy.h:

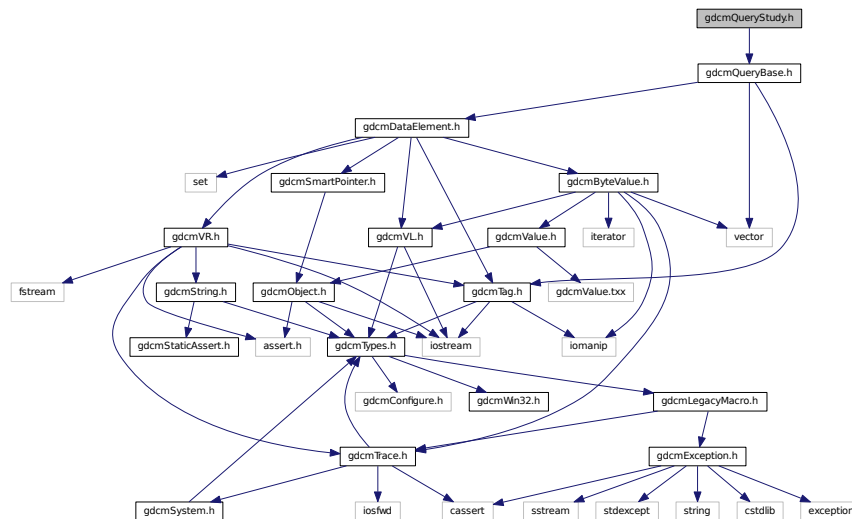

This graph shows which files directly or indirectly include this file:

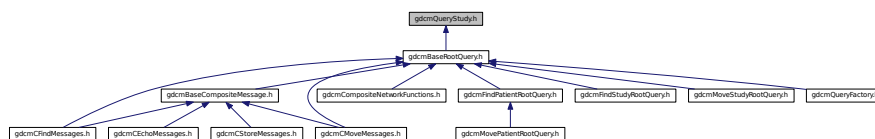

### Classes

- class [gdcm::QueryStudy](#)

*QueryStudy.h contains: class to construct a study-based query for C-FIND and C-MOVE.*


### Namespaces

- [gdcm](#)

## 28.201 gdcmraw.dox File Reference

## 28.202 gdcmRAWCodec.h File Reference

```
#include "gdcmImageCodec.h"
```

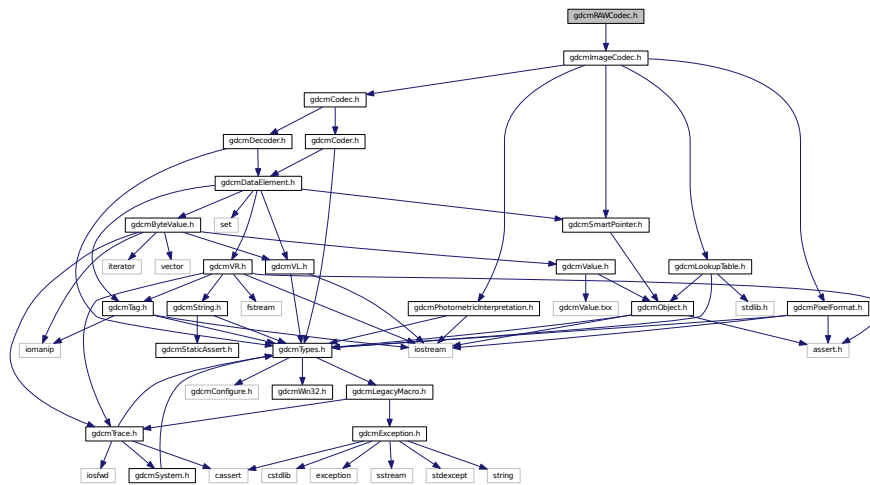

## Classes

- class `gdcm::RAWCodec`

*RAWCodec* class.

## Namespaces

- **gdcm**

## 28.203 gdcmReader.h File Reference

```
#include "gdcmFile.h"
#include <fstream>
```



```

graph BT
    gdcImageRegionReader[gdcImageRegionReader.h] --> gdcImageReader[gdcImageReader.h]
    gdcImageReader --> gdcPixmapReader[gdcPixmapReader.h]
    gdcSurfaceReader[gdcSurfaceReader.h] --> gdcSegmentReader[gdcSegmentReader.h]
    gdcImageRegionReader --> gdcReader[gdcReader.h]
    gdcImageReader --> gdcReader
    gdcPixmapReader --> gdcReader
    gdcSurfaceReader --> gdcReader
    gdcSegmentReader --> gdcReader
    gdcStreamImageReader[gdcStreamImageReader.h] --> gdcReader


```

- class `gdcm::Reader`  
*Reader* ala DOM (Document *Object* Model)


- **gdcm**

## 28.204 gdcmRegion.h File Reference

```
#include "gdcmTypes.h"
#include <vector>
#include <iostream>
Include dependency graph for gdcmRegion.h:
```

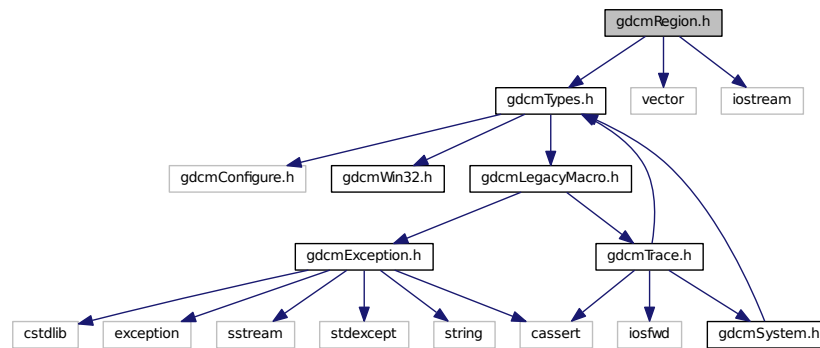

This graph shows which files directly or indirectly include this file:

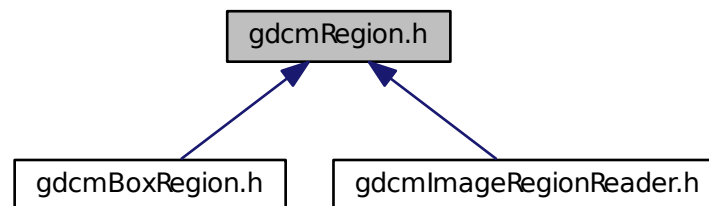

### Classes

- class `gdcm::Region`  
*Class for manipulation region.*

### Namespaces

- `gdcm`

## Functions


- `std::ostream & gdcm::operator<<` (`std::ostream &os`, `const Region &r`)

## 28.205 gdcmRescaler.h File Reference

```
#include "gdcmTypes.h"
```

```
#include "gdcmPixelFormat.h"
```

Include dependency graph for `gdcmRescaler.h`:

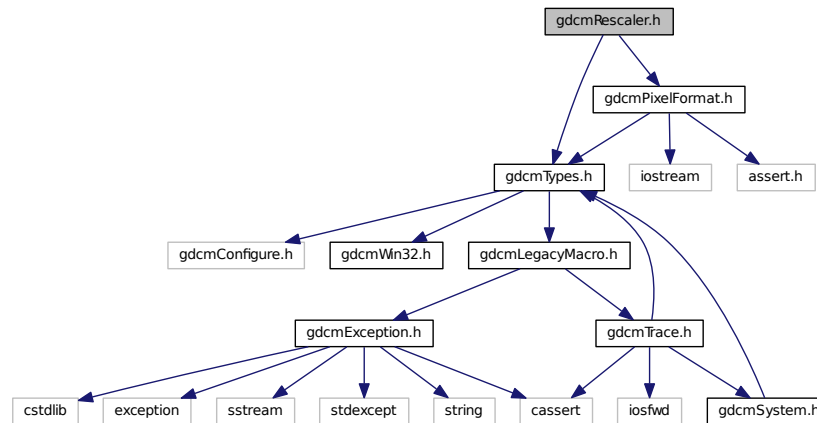

## Classes

- class `gdcm::Rescaler`

*Rescale class* This class is meant to apply the linear transform of Stored Pixel *Value* to Real World *Value*. This is mostly found in CT or PET dataset, where the value are stored using one type, but need to be converted to another scale using a linear transform. There are basically two cases: In CT: the linear transform is generally integer based. E.g. the Stored Pixel *Type* is unsigned short 12bits, but to get Hounsfield unit, one need to apply the linear transform:

$$RWV = 1. * SV - 1024$$

So the best scalar to store the Real World *Value* will be 16 bits signed type.

## Namespaces

- `gdcm`

## 28.206 gdcmRLECodec.h File Reference


```
#include "gdcmImageCodec.h"
```



- Class to do RLE.*

- **gdcm**

```
#include "gdcmTypes.h"
```

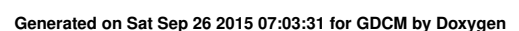

## Classes

- class [gdcm::network::RoleSelectionSub](#)

[RoleSelectionSub](#) PS 3.7 [Table D.3-9 SCP/SCU ROLE SELECTION SUB-ITEM FIELDS \(A-ASSOCIATE-RQ\)](#)

## Namespaces


- [gdcm](#)
- [gdcm::network](#)

## 28.208 gdcmscanner.dox File Reference

## 28.209 gdcmScanner.h File Reference

```
#include "gdcmDirectory.h"
#include "gdcmSubject.h"
#include "gdcmTag.h"
#include "gdcmPrivateTag.h"
#include "gdcmSmartPointer.h"
#include <map>
#include <set>
#include <string>
#include <string.h>
```

Include dependency graph for gdcmScanner.h:

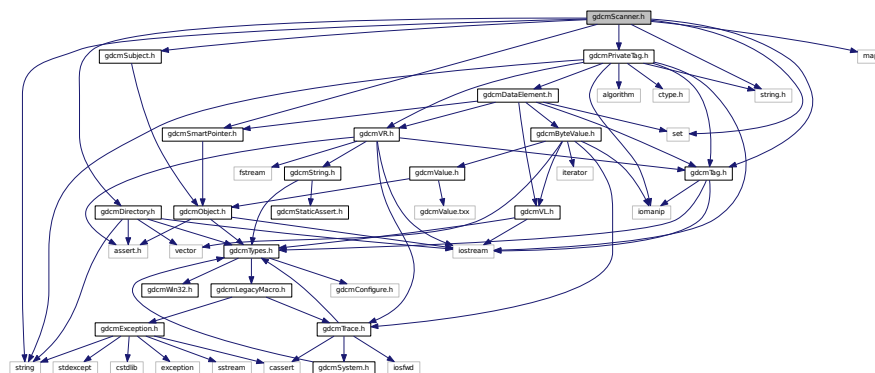

## Classes

- struct [gdcm::Scanner::Itstr](#)
- class [gdcm::Scanner](#)

[Scanner](#) This filter is meant for quickly browsing a [FileSet](#) (a set of files on disk). Special consideration are taken so as to read the minimum amount of information in each file in order to retrieve the user specified set of DICOM [Attribute](#).

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Scanner &s)`


## 28.210 gdcmscu.dox File Reference

## 28.211 gdcmSegment.h File Reference

```
#include <vector>
#include <gdcmObject.h>
#include <gdcmSurface.h>
#include "gdcmSegmentHelper.h"
Include dependency graph for gdcmSegment.h:
```

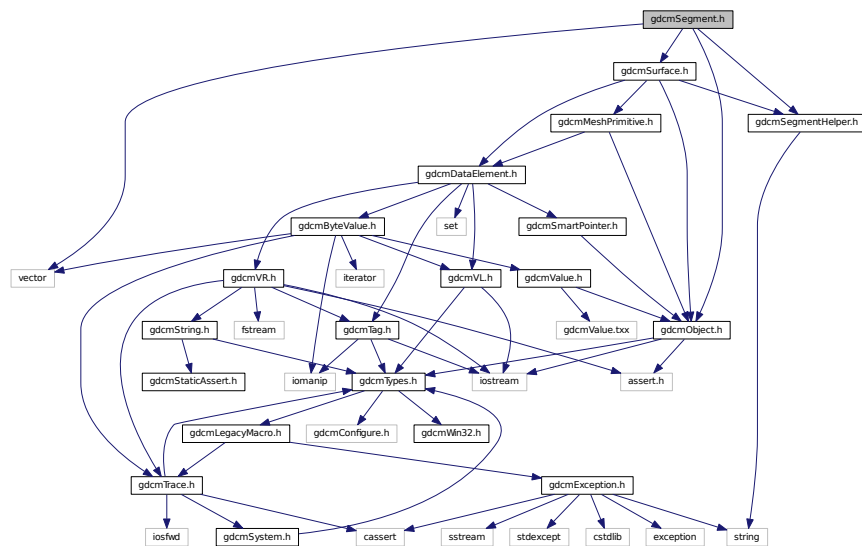

This graph shows which files directly or indirectly include this file:

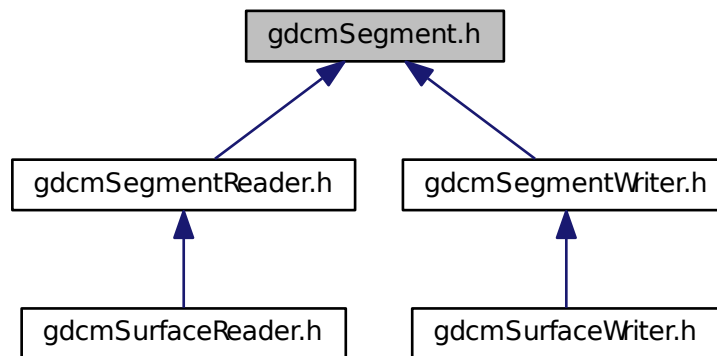

## Classes

- class [gdcm::Segment](#)

*This class defines a segment. It mainly contains attributes of group 0x0062. In addition, it can be associated with surface.*


## Namespaces

- [gdcm](#)

## 28.212 gdcmSegmentedPaletteColorLookupTable.h File Reference

```
#include "gdcmLookupTable.h"
```

Include dependency graph for `gdcmsSegmentedPaletteColorLookupTable.h`:

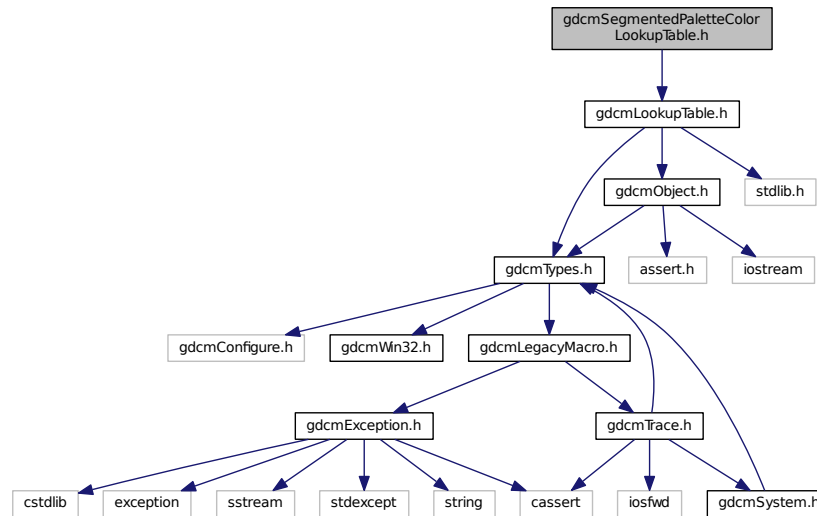

## Classes

- class [gdcms::SegmentedPaletteColorLookupTable](#)  
*SegmentedPaletteColorLookupTable* class.


## Namespaces

- [gdcms](#)

## 28.213 gdcmsSegmentHelper.h File Reference


```
#include <string>
```

Include dependency graph for `gdcmsSegmentHelper.h`:

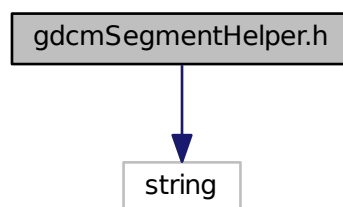



This graph shows which files directly or indirectly include this file:

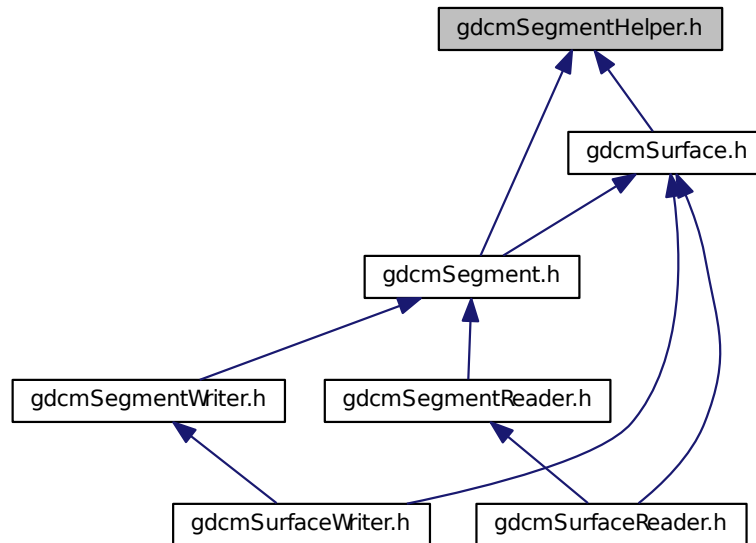

## Classes

- struct [gdcm::SegmentHelper::BasicCodedEntry](#)

*This structure defines a basic coded entry with all of its attributes.*

## Namespaces

- [gdcm](#)
- [gdcm::SegmentHelper](#)


## 28.214 gdcmSegmentReader.h File Reference

```


#include <map>
#include <gdcmReader.h>
#include <gdcmSegment.h>

```

Include dependency graph for `gdcmSegmentReader.h`:

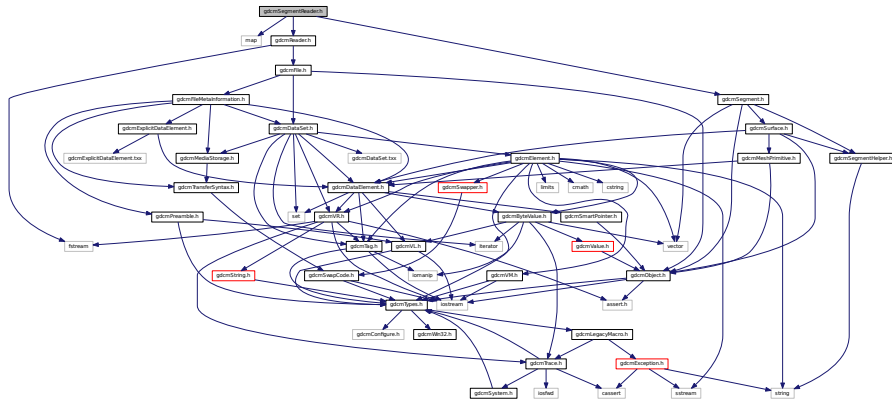

This graph shows which files directly or indirectly include this file:

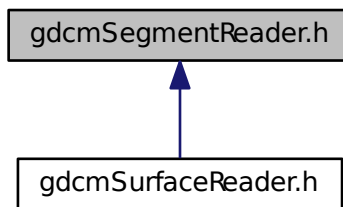

## Classes

- class `gdcm::SegmentReader`

*This class defines a segment reader. It reads attributes of group 0x0062.*

## Namespaces

- `gdcm`

## 28.215 gdcmSegmentWriter.h File Reference

```
#include <gdcmWriter.h>
#include <gdcmSegment.h>
```

```
graph BT; gdcmsurfacewriter[hgdcmsurfaceWriter.h] --> gdcmsegmentwriter[hgdcmsegmentWriter.h];
```


- class `gdcm::SegmentWriter`

## Namespaces

- **gdcm**

```
#include "gdcmValue.h"
```

```
#include "gdcmlVL.h"
#include "gdcmlFragment.h"
#include "gdcmlBasicOffsetTable.h"
Include dependency graph for gdcmlSequenceOfFragments.h:
```

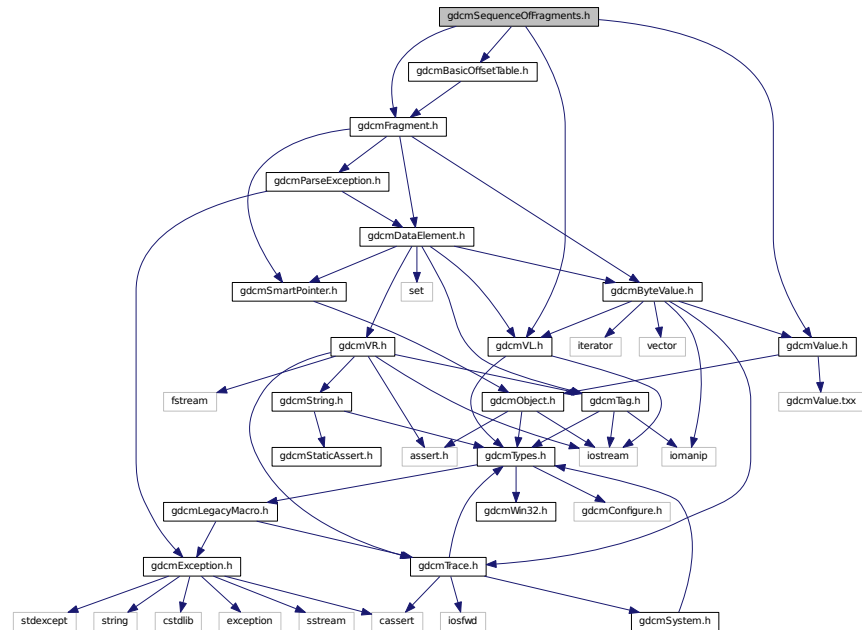

## Classes

- class `gdcmm::SequenceOfFragments`  
*Class to represent a Sequence Of Fragments.*

## Namespaces

- gdc

## 28.217 gdcmSequenceOfItems.h File Reference


```
#include "gdcmValue.h"
#include "gdcmItem.h"
#include <vector>
#include <cstring>
#include "gdcmSequenceOfItems.hxx"
```

- class `gdcmm::SequenceOfItems`  
*Class to represent a Sequence Of Items (value representation : SQ)*

- **gdcm**

```
#include "gdcmTag.h"
#include "gdcmSmartPointer.h"
#include "gdcmFile.h"
#include <vector>
#include <string>
#include <map>
```

**Include dependency graph for gdcmSerieHelper.h:**

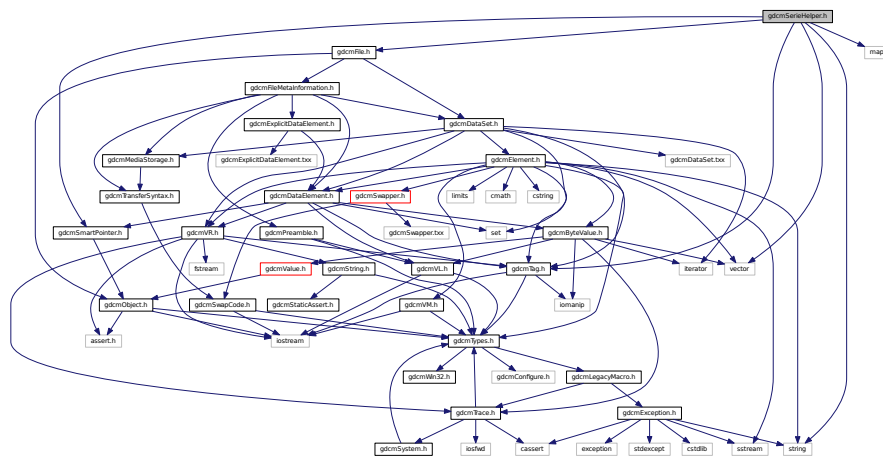

## Classes

- class [gdcmm::FileWithName](#)

*FileWithName.*

- struct [gdcmm::SerieHelper::Rule](#)
- class [gdcmm::SerieHelper](#)

*SerieHelper DO NOT USE this class, it is only a temporary solution for ITK migration from GDCM 1.x to GDCM 2.x It will disappear soon, you've been warned.*

## Namespaces

- [gdcmm](#)

## Typedefs

- typedef bool(\* [gdcmm::BOOL\\_FUNCTION\\_PFILE\\_PFILE\\_POINTER](#)) (File \*, File \*)
- typedef std::vector< SmartPointer< FileWithName > > [gdcmm::FileList](#)


## Enumerations

- enum [gdcmm::CompOperators](#) {  
[gdcmm::GDCM\\_EQUAL](#) = 0,  
[gdcmm::GDCM\\_DIFFERENT](#),  
[gdcmm::GDCM\\_GREATER](#),  
[gdcmm::GDCM\\_GREATEROREQUAL](#),  
[gdcmm::GDCM\\_LESS](#),  
[gdcmm::GDCM\\_LESSEQUAL](#) }
- enum [gdcmm::LodModeType](#) {  
[gdcmm::LD\\_ALL](#) = 0x00000000,  
[gdcmm::LD\\_NOSEQ](#) = 0x00000001,  
[gdcmm::LD\\_NOSHADOW](#) = 0x00000002,  
[gdcmm::LD\\_NOSHADOWSEQ](#) = 0x00000004 }


## 28.219 gdcmmSeries.h File Reference

```
#include "gdcmmTypes.h"
```

Include dependency graph for gdcmSeries.h:

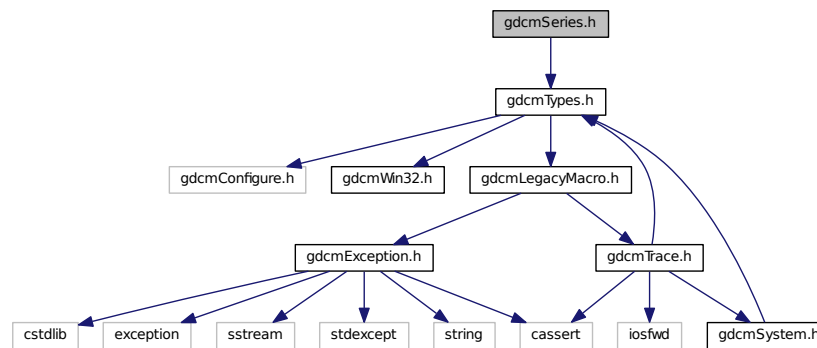

This graph shows which files directly or indirectly include this file:

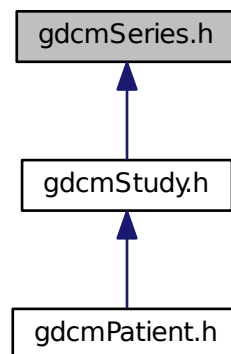

## Classes

- class `gdcm::Series`  
*Series.*


## Namespaces

- `gdcm`


## 28.220 gdcmServiceClassApplicationInformation.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmServiceClassApplicationInformation.h:

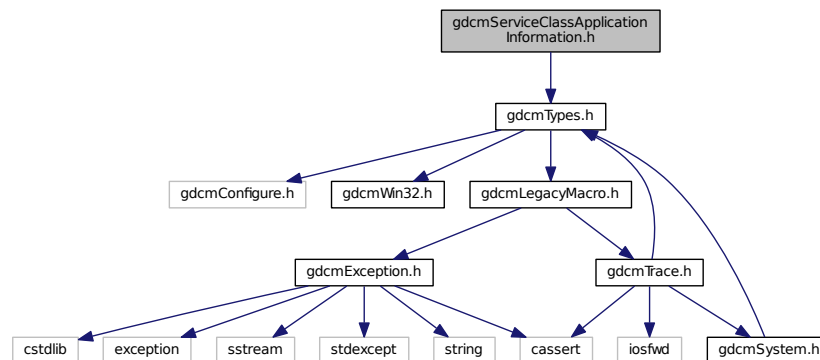

This graph shows which files directly or indirectly include this file:

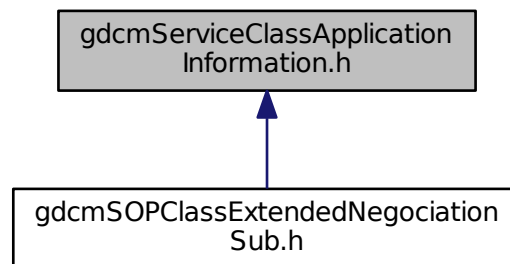

### Classes


- class [gdcm::network::ServiceClassApplicationInformation](#)

### Namespaces

- [gdcm](#)
- [gdcm::network](#)



```
#include "gdcmSubject.h"
#include "gdcmPresentationContext.h"
#include "gdcmFile.h"
#include "gdcmNetworkStateID.h"
Include dependency graph for gdcmServiceClassUser.h:
```

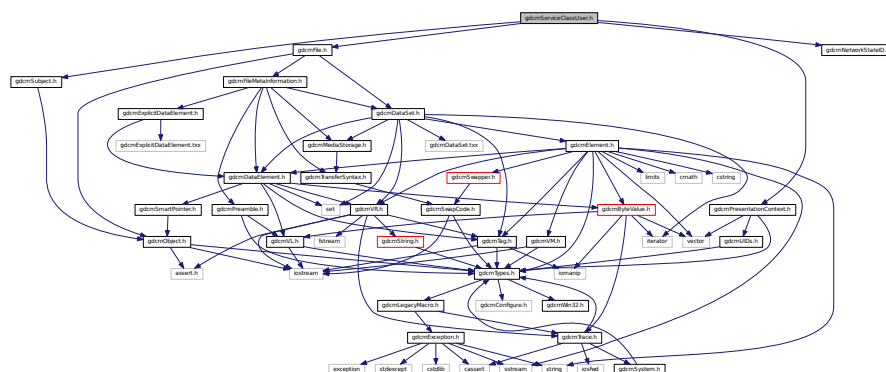


- class `gdcm::ServiceClassUser`

*ServiceClassUser.*

- `gdcm`
- `gdcm::network`

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmSHA1.h`:

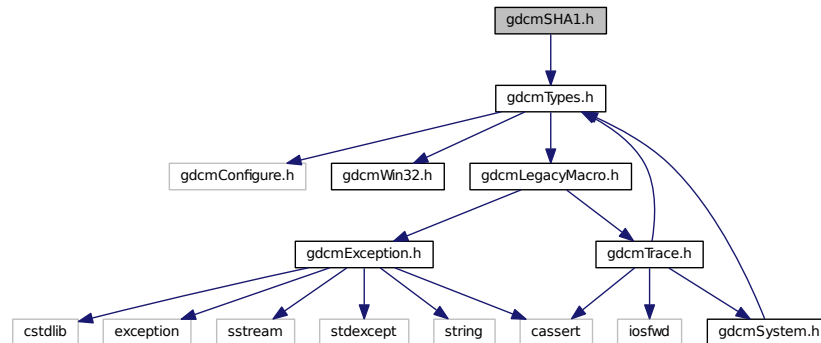

## Classes

- class `gdcm::SHA1`

*Class for `SHA1`.*

## Namespaces

- `gdcm`


## 28.223 `gdcmSimpleSubjectWatcher.h` File Reference

```

#include "gdcmSubject.h"
#include "gdcmCommand.h"
#include "gdcmSmartPointer.h"
#include "gdcmAnonymizeEvent.h"
#include "gdcmDataEvent.h"

```

Include dependency graph for gdcmSimpleSubjectWatcher.h:

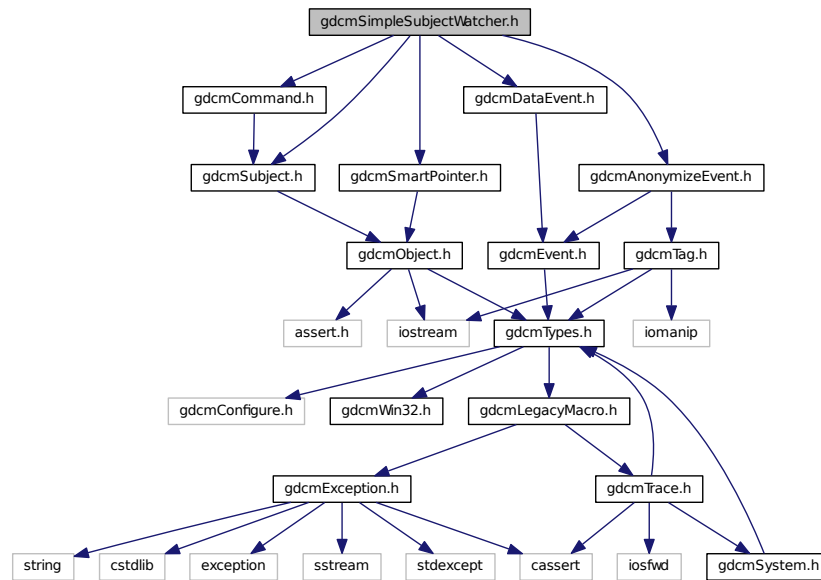

## Classes

- class [gdcm::SimpleSubjectWatcher](#)



*[SimpleSubjectWatcher](#) This is a typical [Subject](#) Watcher class. It will observe all events.*

## Namespaces

- [gdcm](#)

## 28.224 gdcmSmartPointer.h File Reference

```
#include "gdcmObject.h"
```

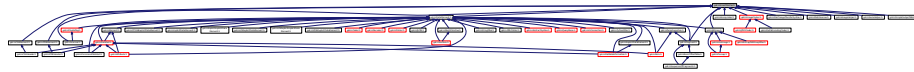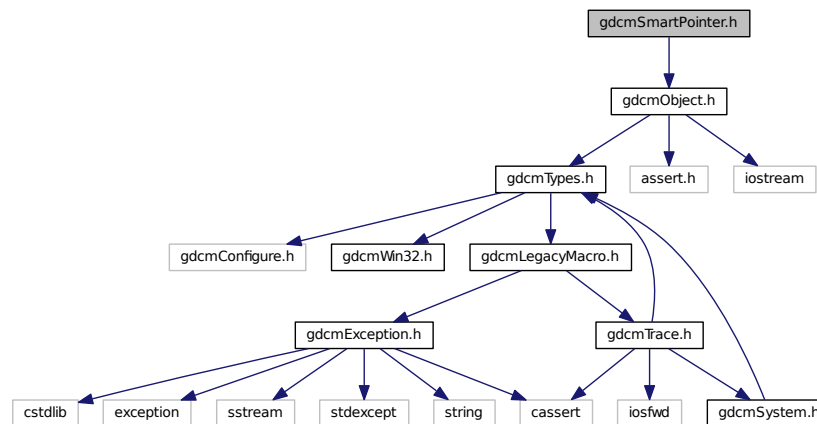


- class `gdcm::SmartPointer< ObjectType >`

- **gdcm**

```
#include "gdcmServiceClassApplicationInformation.h"
```

```
#include "gdcmServiceClassApplicationInformation.h"
```

Include dependency graph for gdcmSOPClassExtendedNegociationSub.h:

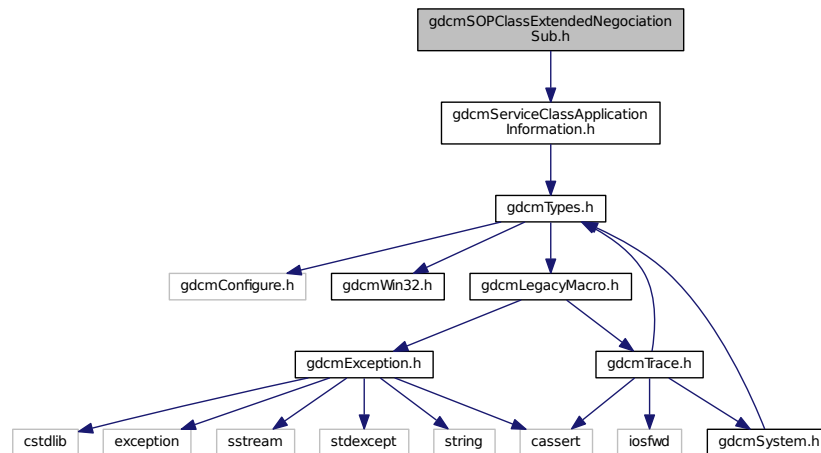

## Classes

- class [gdcm::network::SOPClassExtendedNegociationSub](#)

[SOPClassExtendedNegociationSub](#) PS 3.7 Table D.3-11 SOP CLASS EXTENDED NEGOTIATION SUB-ITEM FIELDS (A-ASSOCIATE-RQ and A-ASSOCIATE-AC)


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.226 gdcmSOPClassUIDToIOD.h File Reference

```
#include "gdcmUIDs.h"
```

Include dependency graph for `gdcmSOPClassUIDToIOD.h`:

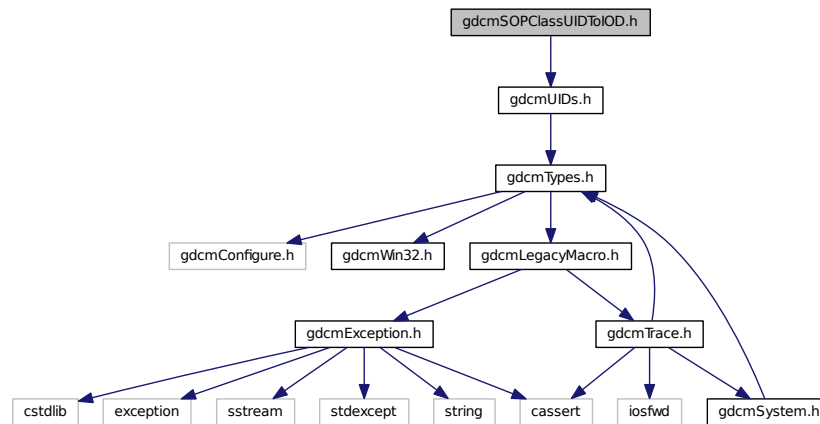

## Classes

- class `gdcm::SOPClassUIDToIOD`

*Class convert a class SOP Class UID into [IOD](#).*

## Namespaces

- `gdcm`


## 28.227 gdcmSorter.h File Reference

```


#include "gdcmDirectory.h"
#include "gdcmTag.h"
#include <vector>
#include <string>
#include <map>

```

Include dependency graph for gdcmSorter.h:

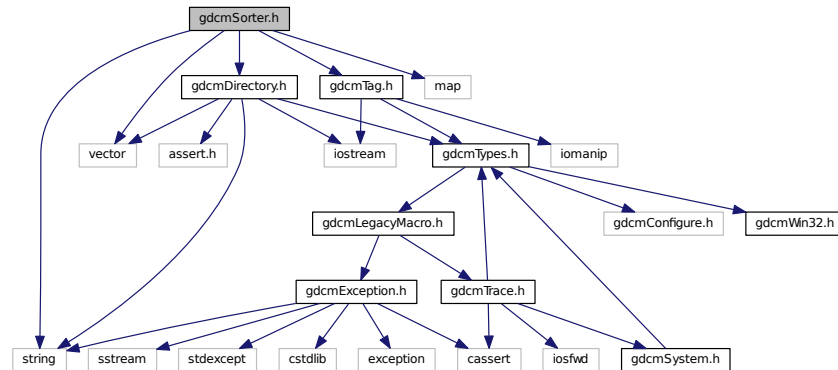

This graph shows which files directly or indirectly include this file:

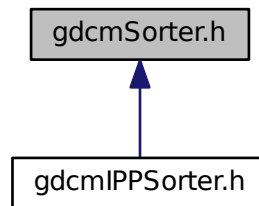

## Classes

- class [gdcm::Sorter](#)

[Sorter](#) General class to do sorting using a custom function You simply need to provide a function of type: [Sorter::Sort↔Function](#).

## Namespaces

- [gdcm](#)


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &os, const Sorter &s)`





Include dependency graph for gdcmSpectroscopy.h:

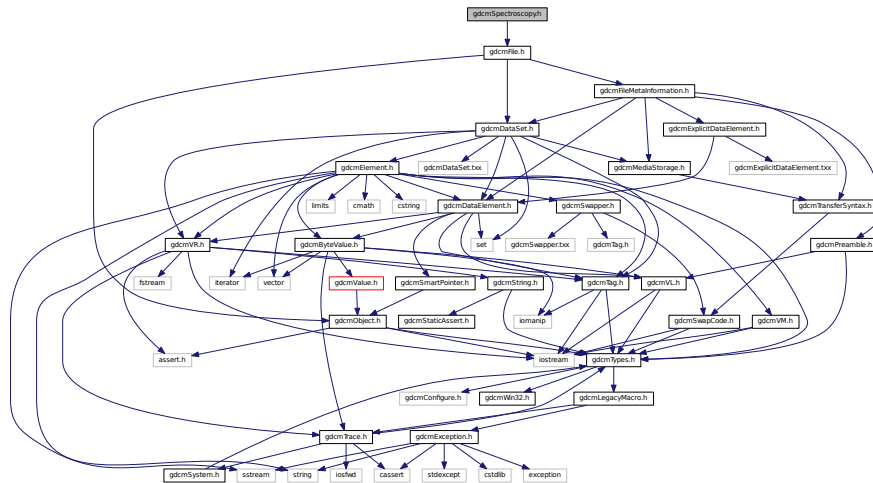

## Classes

- class `gdcm::Spectroscopy`  
*Spectroscopy* class.

## Namespaces


- `gdcm`

## 28.230 gdcmSplitMosaicFilter.h File Reference

```
#include "gdcmFile.h"
```

```
#include "gdcmImage.h"
```

Include dependency graph for `gdcmSplitMosaicFilter.h`:

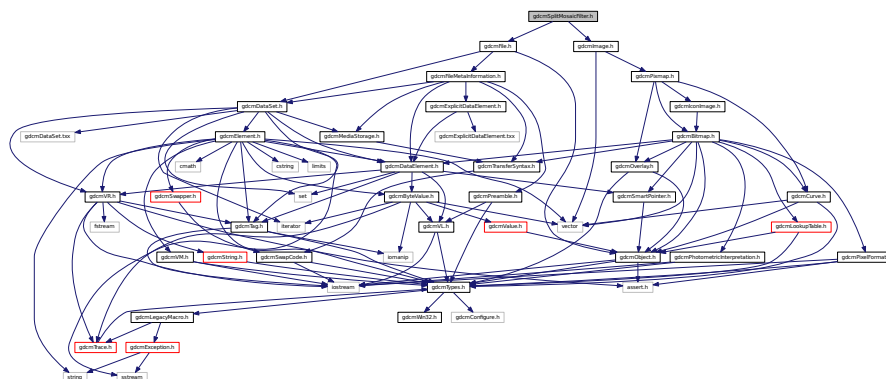

## Classes

- class [gdcm::SplitMosaicFilter](#)


[SplitMosaicFilter](#) class Class to reshuffle bytes for a SIEMENS Mosaic image Siemens CSA [Image](#) Header CSA:= Common Siemens Architecture, sometimes also known as Common syngo Architecture.

## Namespaces

- [gdcm](#)

## 28.231 gdcmStaticAssert.h File Reference

This graph shows which files directly or indirectly include this file:

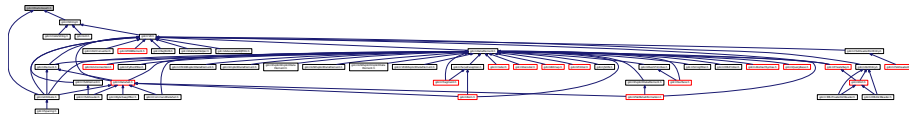

## Classes

- struct [gdcm::static\\_assert\\_test< x >](#)
- struct [gdcm::STATIC\\_ASSERTION\\_FAILURE< x >](#)
- struct [gdcm::STATIC\\_ASSERTION\\_FAILURE< true >](#)

## Namespaces

- [gdcm](#)

## Macros

- #define [GDCM\\_DO\\_JOIN\(X, Y\) GDCM\\_DO\\_JOIN2\(X,Y\)](#)
- #define [GDCM\\_DO\\_JOIN2\(X, Y\) X##Y](#)
- #define [GDCM\\_JOIN\(X, Y\) GDCM\\_DO\\_JOIN\( X, Y \)](#)
- #define [GDCM\\_STATIC\\_ASSERT\(B\)](#)

The [GDCM\\_JOIN](#) + **LINE** is needed to create a uniq identifier.

### 28.231.1 Macro Definition Documentation

28.231.1.1 #define [GDCM\\_DO\\_JOIN\( X, Y \) GDCM\\_DO\\_JOIN2\(X,Y\)](#)


28.231.1.2 #define [GDCM\\_DO\\_JOIN2\( X, Y \) X##Y](#)

28.231.1.3 #define [GDCM\\_JOIN\( X, Y \) GDCM\\_DO\\_JOIN\( X, Y \)](#)

**Value:**

The GDCM\_JOIN + **LINE** is needed to create a uniq identifier.

```
#include "gdcmReader.h"
Include dependency graph for gdcmStreamImageReader.h:
```

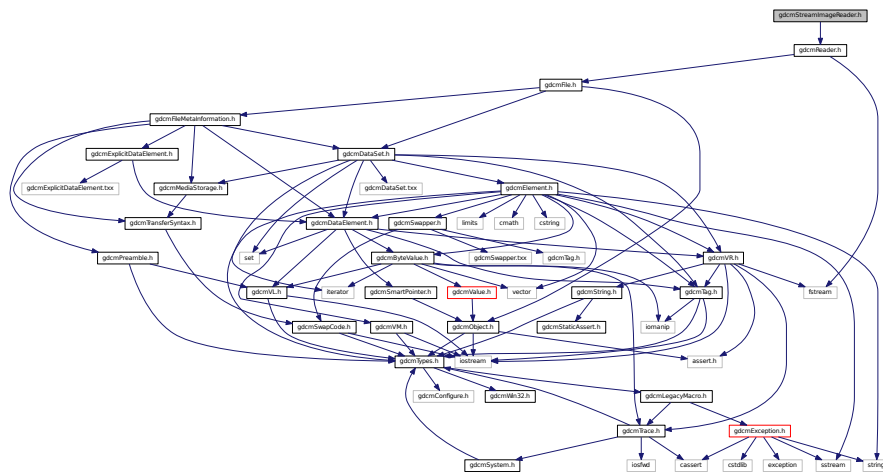

- class `gdcm::StreamImageReader`  
*StreamImageReader*.


- **gdcm**

```
#include "gdcmWriter.h"
#include <iostream>
#include "gdcmDataSet.h"
```


- class `gdcm::StreamImageWriter`  
*StreamImageReader*.

- **gdcm**

```
#include "gdcmTypes.h"
#include "gdcmStaticAssert.h"
Include dependency graph for gdcmString.h:
```

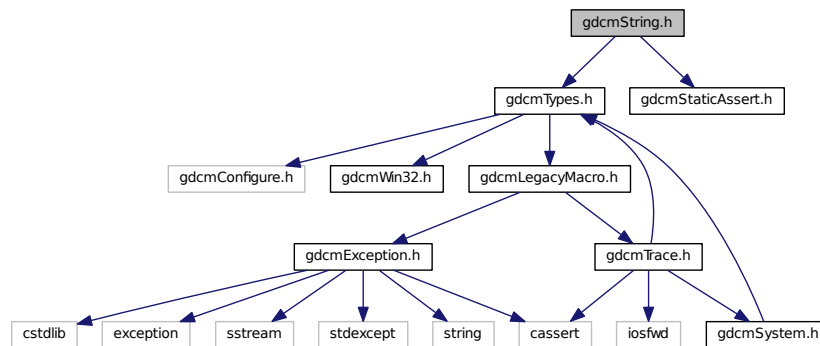

This graph shows which files directly or indirectly include this file:

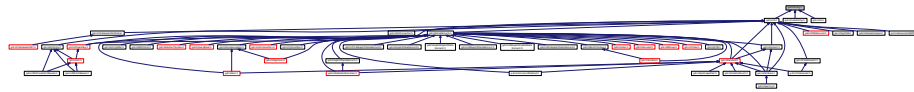

## Classes

- class [gdcm::String< TDelimiter, TMaxLength, TPadChar >](#)  
*String.*

## Namespaces

- [gdcm](#)


## Functions

- template<char TDelimiter, unsigned int TMaxLength, char TPadChar>  
std::istream & [gdcm::operator>>](#) (std::istream &is, String< TDelimiter, TMaxLength, TPadChar > &ms)

## 28.235 gdcmStringFilter.h File Reference

```
#include "gdcmDataElement.h"
#include "gdcmDicts.h"
#include "gdcmFile.h"
```

Include dependency graph for gdcmStringFilter.h:

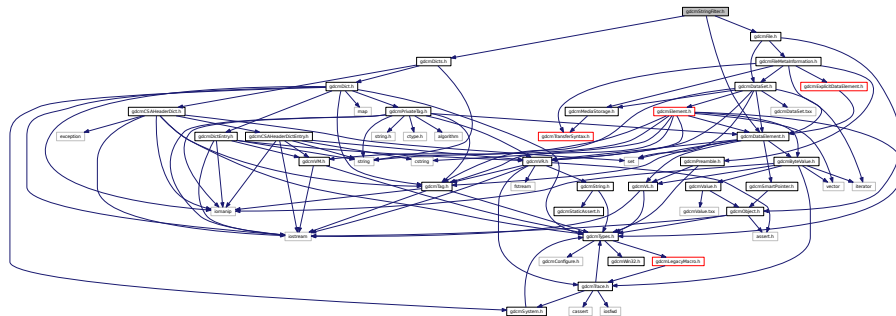

## Classes

- class [gdcm::StringFilter](#)  
*StringFilter* [StringFilter](#) is the class that make *gdcm2.x* looks more like *gdcm1* and transform the binary blob contained in a [DataElement](#) into a string, typically this is a nice feature to have for wrapped language.


## Namespaces

- [gdcm](#)


## 28.236 gdcmStudy.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmSeries.h"
#include <vector>
```

Include dependency graph for gdcmStudy.h:

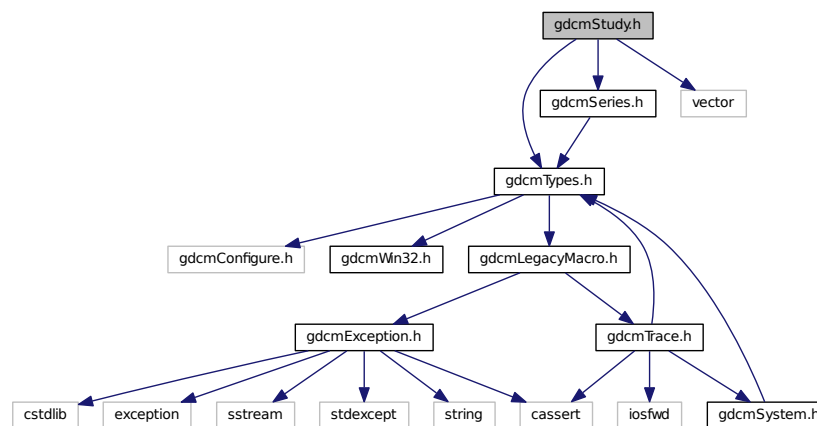

This graph shows which files directly or indirectly include this file:

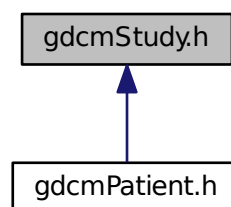

## Classes

- class [gdcm::Study](#)  
*Study.*


## Namespaces

- [gdcm](#)


## 28.237 gdcmSubject.h File Reference

```
#include "gdcmObject.h"
```

Include dependency graph for gdcmSubject.h:

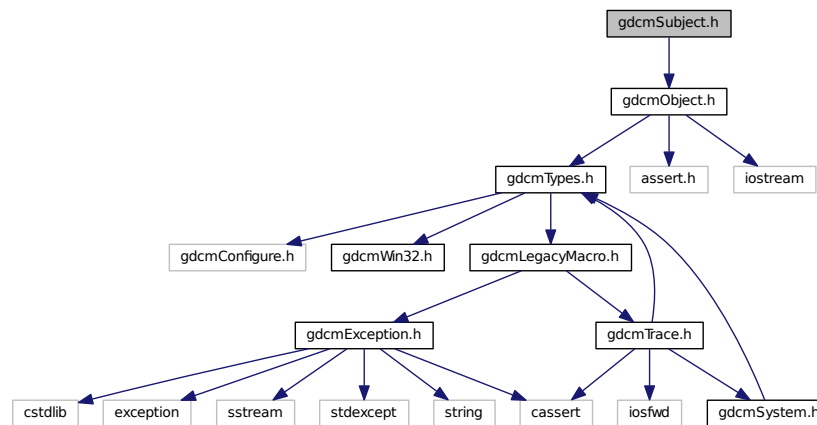

This graph shows which files directly or indirectly include this file:

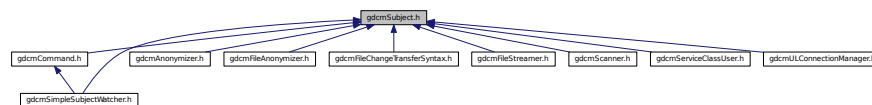

## Classes

- class [gdcm::Subject](#)  
*Subject.*

## Namespaces

- [gdcm](#)

## 28.238 gdcmSurface.h File Reference

```
#include <gdcmObject.h>
```






## 28.239 gdcmSurfaceHelper.h File Reference

```
#include "gdcmTypes.h"
```

```
#include <vector>
```

```
#include <iostream>
```

Include dependency graph for gdcmSurfaceHelper.h:

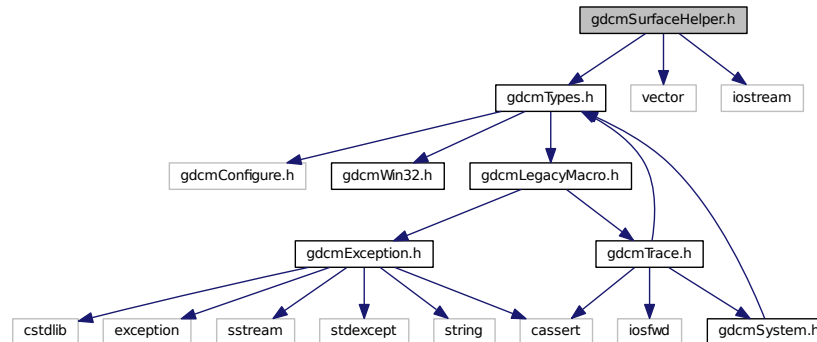

## Classes

- class [gdcm::SurfaceHelper](#)

*[SurfaceHelper](#) Helper class for [Surface](#) object.*


## Namespaces

- [gdcm](#)

## 28.240 gdcmSurfaceReader.h File Reference

```
#include <gdcmSegmentReader.h>
```

```
#include <gdcmSurface.h>
```

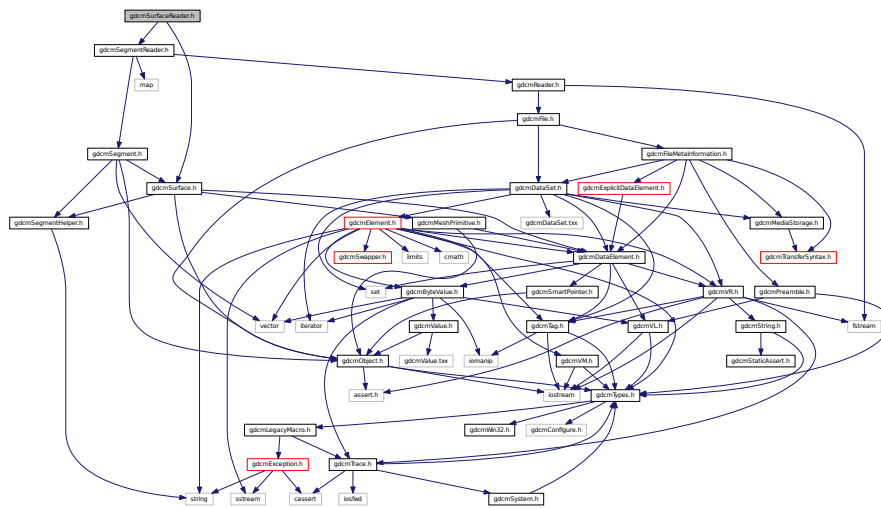

- class `gdcm::SurfaceReader`


*This class defines a SURFACE IE reader. It reads surface mesh module attributes.*

- **gdcm**

## 28.241 gdcmSurfaceWriter.h File Reference

```
#include <gdcmSegmentWriter.h>
#include <gdcmSurface.h>
```

Include dependency graph for gdcmSurfaceWriter.h:

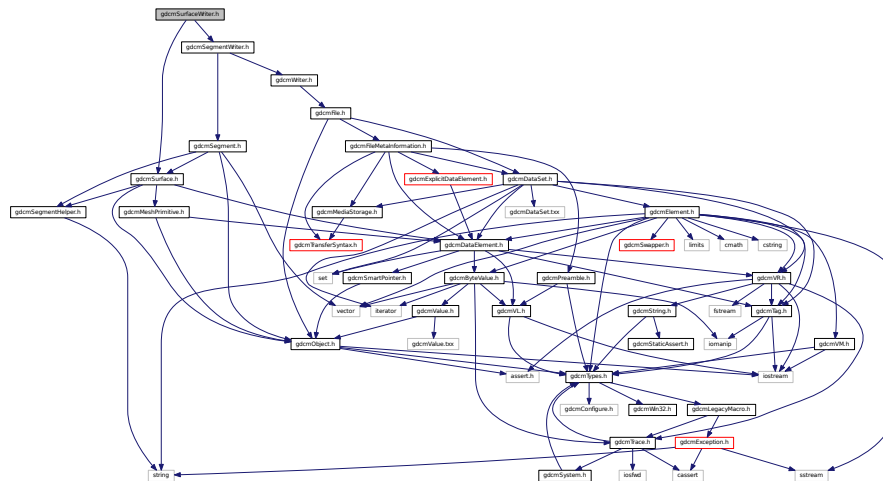

## Classes

- class [gdcm::SurfaceWriter](#)

*This class defines a SURFACE IE writer. It writes surface mesh module attributes.*

## Namespaces


- [gdcm](#)

## 28.242 gdcmSwapCode.h File Reference


```
#include "gdcmTypes.h"
```

```
#include <iostream>
```

Include dependency graph for gdcmSwapCode.h:

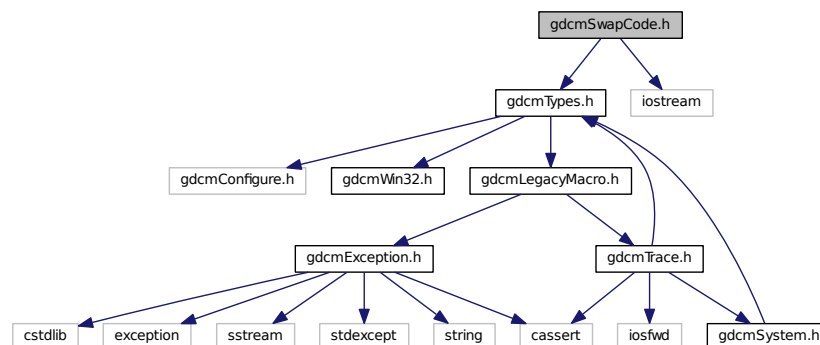

This graph shows which files directly or indirectly include this file:

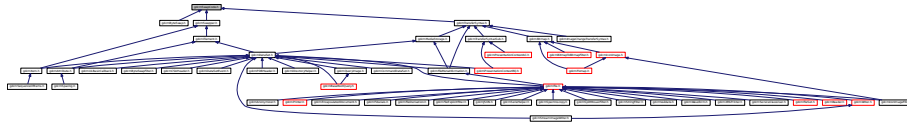

## Classes

- class [gdcm::SwapCode](#)  
*SwapCode* representation.

## Namespaces

- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const SwapCode &sc)`

## 28.243 gdcmSwapper.h File Reference


```
#include "gdcmSwapCode.h"
```

```
#include "gdcmSwapper.txx"
```

Include dependency graph for `gdcmSwapper.h`:

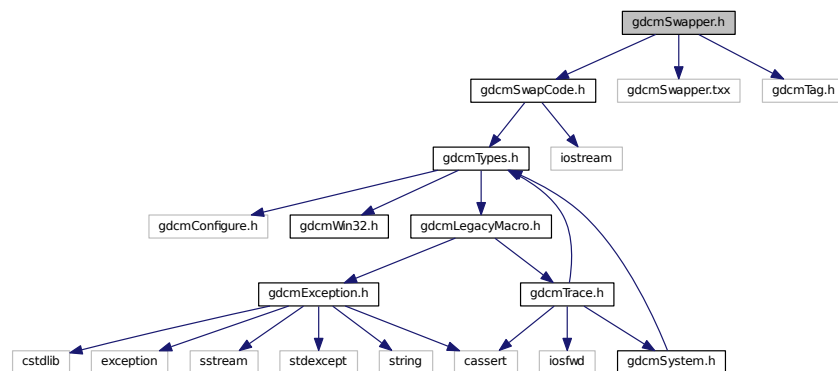

This graph shows which files directly or indirectly include this file:

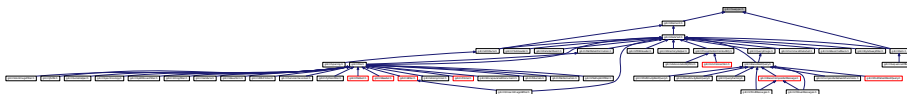

## Classes

- class [gdcm::SwapperDoOp](#)
- class [gdcm::SwapperNoOp](#)


## Namespaces

- [gdcm](#)


## 28.244 gdcmSystem.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmSystem.h:

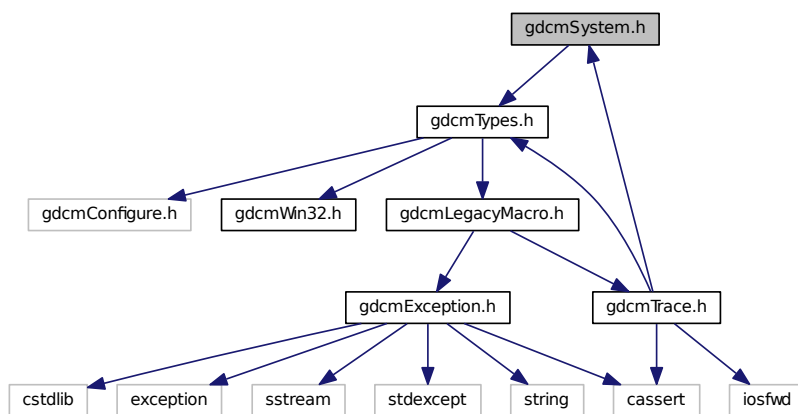

This graph shows which files directly or indirectly include this file:

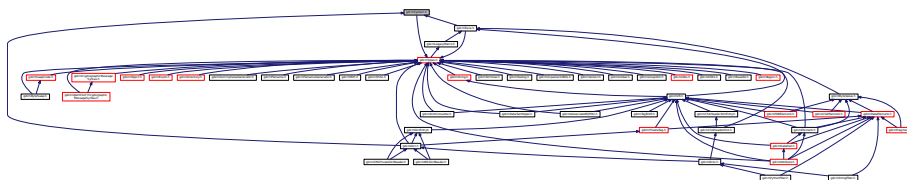

## Classes

- class [gdcm::System](#)  
*Class to do system operation.*


## Namespaces

- [gdcm](#)

## 28.245 gdcmTable.h File Reference

```
#include "gdcmTableEntry.h"
#include "gdcmTag.h"
#include <map>
```

Include dependency graph for gdcmTable.h:

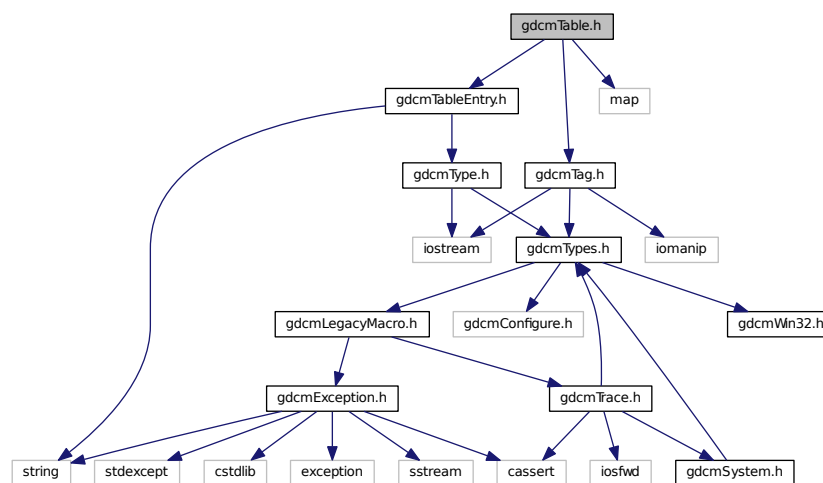

## Classes

- class [gdcm::Table](#)  
*Table.*


## Namespaces

- [gdcm](#)


## 28.246 gdcmTableEntry.h File Reference

```
#include "gdcmType.h"
#include <string>
```

Include dependency graph for gdcmTableEntry.h:

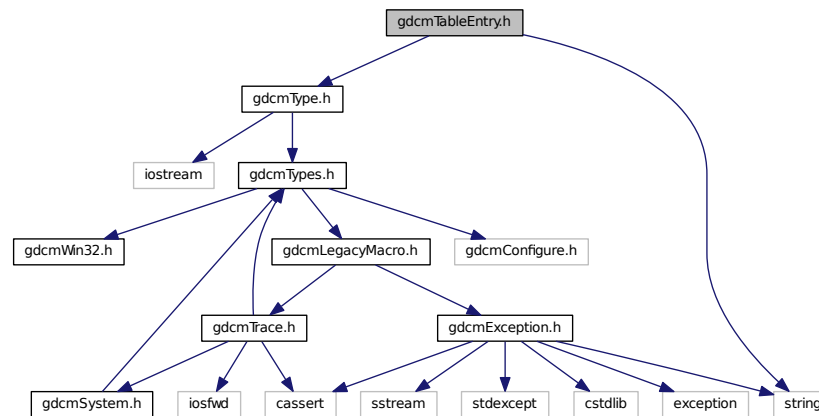

This graph shows which files directly or indirectly include this file:

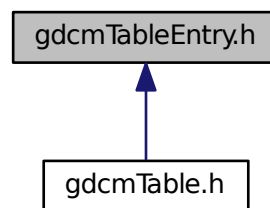

## Classes

- class `gdcm::TableEntry`  
*TableEntry*.


## Namespaces

- `gdcm`

## 28.247 gdcmTableReader.h File Reference


```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmTableReader.h:

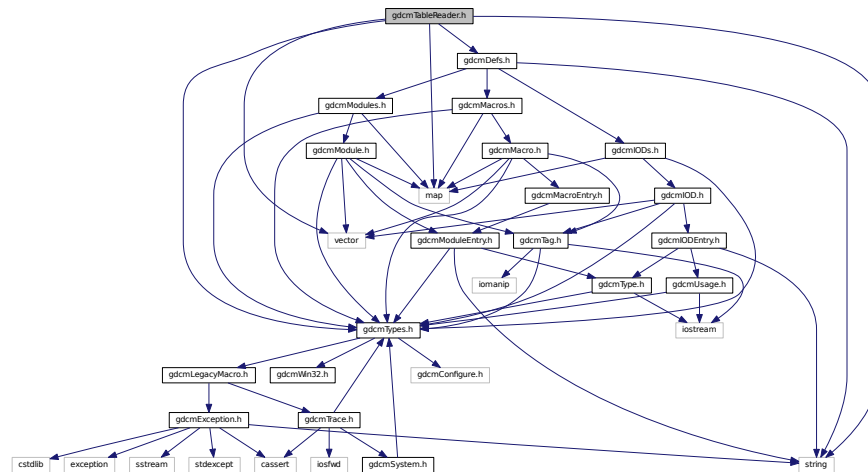

```
graph BT; gdcXMLDictReader.h --> gdcTableReader.h; gdcXMLPrivateDictReader.h --> gdcTableReader.h
```

The diagram illustrates the inheritance structure of the GDC library. At the top is a box labeled `gdcTableReader.h`. Below it are two boxes: `gdcXMLDictReader.h` on the left and `gdcXMLPrivateDictReader.h` on the right. Blue arrows point from each of these bottom boxes up to the `gdcTableReader.h` box, indicating that both `gdcXMLDictReader.h` and `gdcXMLPrivateDictReader.h` inherit from `gdcTableReader.h`.

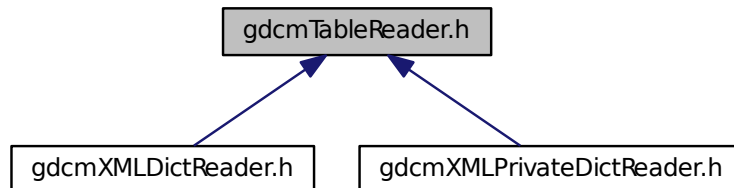


- class `gdcm::TableReader`  
*Class for representing a `TableReader`.*

- [gdc](#)




## 28.248 gdcmTag.h File Reference

```
#include "gdcmTypes.h"
#include <iostream>
#include <iomanip>
Include dependency graph for gdcmTag.h:
```

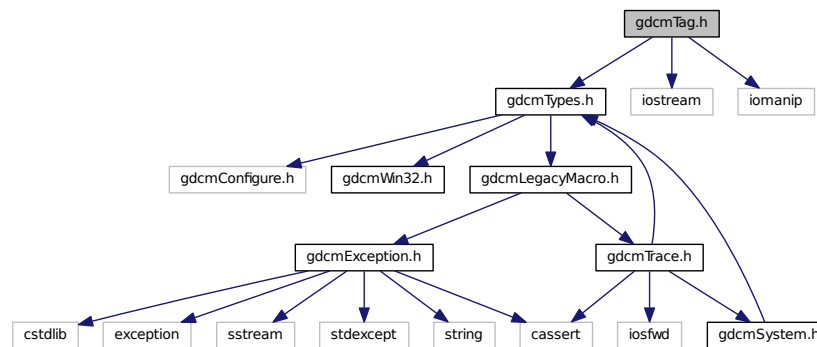

This graph shows which files directly or indirectly include this file:

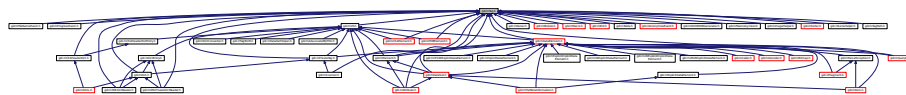

### Classes

- class [gdcm::Tag](#)

Class to represent a DICOM Data [Element](#) ([Attribute](#)) [Tag](#) ([Group](#), [Element](#)). Basically an `uint32_t` which can also be expressed as two `uint16_t` ([group](#) and [element](#))

### Namespaces

- [gdcm](#)

### Functions


- `std::ostream & gdcm::operator<< (std::ostream &_os, const Tag &_val)`
- `std::istream & gdcm::operator>> (std::istream &_is, Tag &_val)`

## 28.249 gdcmTagPath.h File Reference

```
#include "gdcmTag.h"
```

```
#include <vector>
```

Include dependency graph for `gdcmTagPath.h`:

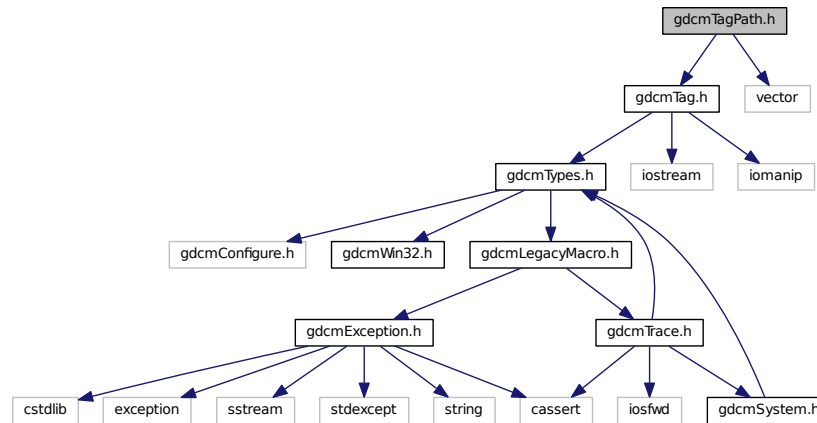

## Classes

- class [gdcm::TagPath](#)

*class to handle a path of tag.*


## Namespaces

- [gdcm](#)

## 28.250 gdcmTagToVR.h File Reference

```
#include "gdcmVR.h"
```

Include dependency graph for gdcmTagToVR.h:

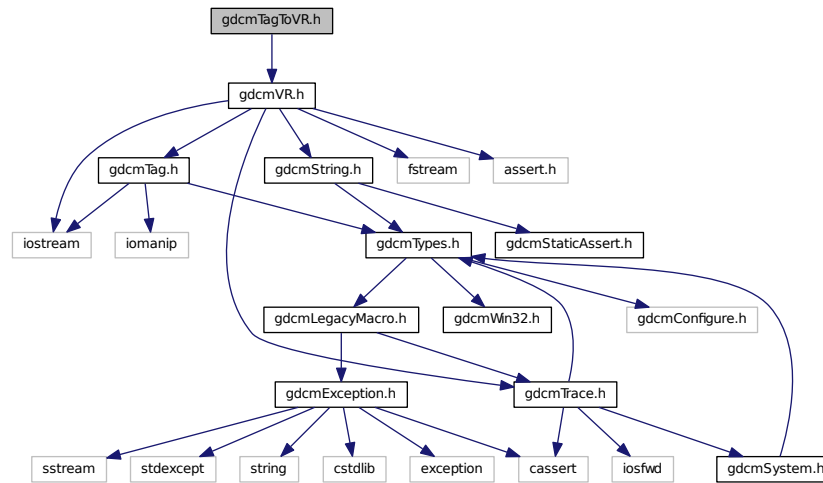

## Namespaces

- [gdcm](#)

## Functions


- VR::VRType [gdcm::GetVRFromTag](#) (Tag const &tag)

## 28.251 gdcmtar.dox File Reference

## 28.252 gdcmTerminal.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcTerminal.h`:

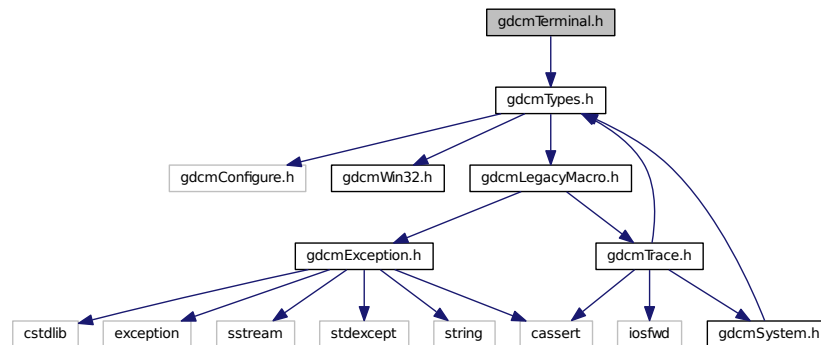

## Namespaces

- `gdc`
- `gdc::terminal`

*Class for Terminal Allow one to print in color in a shell.*

## Enumerations

- enum `gdc::terminal::Attribute` {  
`gdc::terminal::reset` = 0,  
`gdc::terminal::bright` = 1,  
`gdc::terminal::dim` = 2,  
`gdc::terminal::underline` = 3,  
`gdc::terminal::blink` = 5,  
`gdc::terminal::reverse` = 7,  
`gdc::terminal::hidden` = 8 }
- enum `gdc::terminal::Color` {  
`gdc::terminal::black` = 0,  
`gdc::terminal::red`,  
`gdc::terminal::green`,  
`gdc::terminal::yellow`,  
`gdc::terminal::blue`,  
`gdc::terminal::magenta`,  
`gdc::terminal::cyan`,  
`gdc::terminal::white` }
- enum `gdc::terminal::Mode` {  
`gdc::terminal::CONSOLE` = 0,  
`gdc::terminal::VT100` }

## Functions

- `GDCM_EXPORT` `std::string gdc::terminal::setattr` (Attribute att)
- `GDCM_EXPORT` `std::string gdc::terminal::setbgcolor` (Color c)


- [GDCM\\_EXPORT](#) `std::string` [gdcm::terminal::setfgcolor](#) (Color c)
- [GDCM\\_EXPORT](#) `void` [gdcm::terminal::setmode](#) (Mode m)

## 28.253 gdcmTestDriver.h File Reference

```
#include <clocale>
```

```
#include <locale>
```

Include dependency graph for gdcmTestDriver.h:

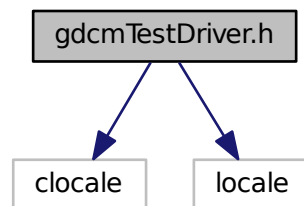


## 28.254 gdcmTesting.h File Reference

```
#include "gdcmTypes.h"
```

```
#include <iostream>
```

Include dependency graph for gdcmTesting.h:

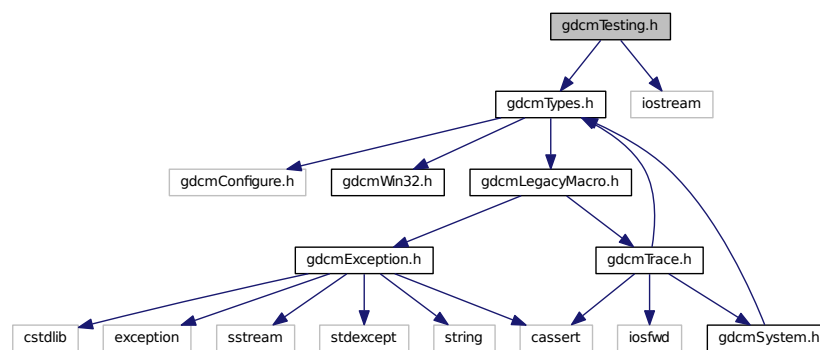

## Classes

- class [gdcm::Testing](#)  
*class for testing*


## Namespaces

- [gdcm](#)


## 28.255 gdcmTrace.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmSystem.h"
#include <iosfwd>
#include <cassert>
```

Include dependency graph for gdcmTrace.h:

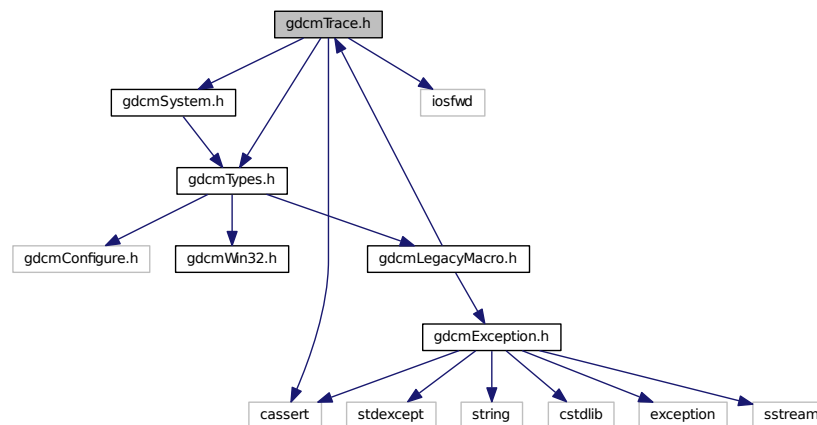

This graph shows which files directly or indirectly include this file:

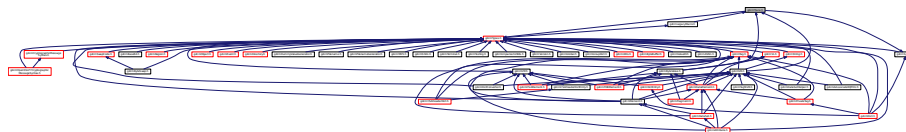

## Classes

- class [gdcm::Trace](#)  
*Trace.*

## Namespaces

- [gdcm](#)

## Macros

- #define `GDCM_FUNCTION` "<unknown>"
- #define `gdcmAssertAlwaysMacro`(arg) `gdcmAssertMacro`(arg)  
*AssertAlways.*
- #define `gdcmAssertMacro`(arg)  
*Assert.*
- #define `gdcmDebugMacro`(msg)  
*Debug.*
- #define `gdcmErrorMacro`(msg)  
*Error this is pretty bad, more than just warning It could mean lost of data, something not handle...*
- #define `gdcmWarningMacro`(msg)  
*Warning.*

### 28.255.1 Macro Definition Documentation

28.255.1.1 #define `GDCM_FUNCTION` "<unknown>"

28.255.1.2 #define `gdcmAssertAlwaysMacro`( *arg* ) `gdcmAssertMacro`(arg)

*AssertAlways.*

Parameters

<i>arg</i>	argument to test An easy solution to pass also a message is to do: <code>gdcmAssertMacro( "my message" &amp;&amp; 2 &lt; 3 )</code>
------------	---

Referenced by `gdcm::SequenceOfFragments::ReadValue()`, and `gdcm::VR::Write()`.

28.255.1.3 #define `gdcmAssertMacro`( *arg* )

Value:

```
{
    if( !(arg) )
    {
        std::ostringstream osmacro;
        osmacro << "Assert: In " __FILE__ ", line " << __LINE__
        << ", function " << GDCM_FUNCTION
        << "\n\n";
        std::ostream &_os = gdcm::Trace::GetErrorStream();
        _os << osmacro.str() << std::endl;
        assert ( arg );
    }
}
```

*Assert.*

Parameters

<i>arg</i>	argument to test An easy solution to pass also a message is to do: <code>gdcmAssertMacro( "my message" &amp;&amp; 2 &lt; 3 )</code>
------------	---

Referenced by `gdcm::PixelFormat::SetSamplesPerPixel()`.

**28.255.1.4 #define gdcDebugMacro( msg )****Value:**

```

{
    if( gdc::Trace::GetDebugFlag() )
    {
        std::ostringstream osmacro;
        osmacro << "Debug: In " __FILE__ ", line " << __LINE__
            << ", function " << GDCM_FUNCTION << '\n'
            << "Last system error was: "
            << gdc::System::GetLastSystemError() << '\n' << msg;
        std::ostream &_os = gdc::Trace::GetDebugStream();
        _os << osmacro.str() << "\n\n" << std::endl;
    }
}

```

Debug.

**Parameters**

<i>msg</i>	message part
------------	--------------

Referenced by `gdc::ByteValue::ByteValue()`, `gdc::OpenSSLCryptoFactory::OpenSSLCryptoFactory()`, `gdc::OpenSSLP7CryptoFactory::OpenSSLP7CryptoFactory()`, `gdc::BasicOffsetTable::Read()`, `gdc::Item::Read()`, `gdc::SequenceOfItems::Read()`, `gdc::VR::Read()`, `gdc::SequenceOfFragments::ReadPreValue()`, `gdc::SequenceOfFragments::ReadValue()`, and `gdc::ByteValue::SetLength()`.

**28.255.1.5 #define gdcErrorMacro( msg )****Value:**

```

{
    if( gdc::Trace::GetErrorFlag() )
    {
        std::ostringstream osmacro;
        osmacro << "Error: In " __FILE__ ", line " << __LINE__
            << ", function " << GDCM_FUNCTION << '\n'
            << msg << "\n\n";
        std::ostream &_os = gdc::Trace::GetErrorStream();
        _os << osmacro.str() << std::endl;
    }
}

```

Error this is pretty bad, more than just warning It could mean lost of data, something not handle...

**Parameters**

<i>msg</i>	second message part
------------	---------------------

Referenced by `gdc::CommandDataSet::Insert()`, `gdc::FileMetaInformation::Insert()`, `gdc::DataSet::Insert()`, `gdc::Item::Read()`, and `gdc::Fragment::ReadBacktrack()`.

**28.255.1.6 #define gdcWarningMacro( msg )****Value:**

```

{
    if( gdc::Trace::GetWarningFlag() )
    {
        std::ostringstream osmacro;
        osmacro << "Warning: In " __FILE__ ", line " << __LINE__
            << ", function " << GDCM_FUNCTION << "\n"
            << msg << "\n\n";
    }
}

```



```

std::ostream &_os = gdcm::Trace::GetWarningStream();
_os << omacro.str() << std::endl;
}
}

```

Warning.

Parameters


<i>msg</i>	message part
------------	--------------

Referenced by `gdcm::DataSet::InsertDataElement()`, `gdcm::Item::Read()`, `gdcm::SequenceOfItems::Read()`, `gdcm::Fragment::ReadBacktrack()`, `gdcm::Fragment::ReadValue()`, `gdcm::SequenceOfFragments::ReadValue()`, `gdcm::OpenSSL7CryptographicMessageSyntax::SetPassword()`, and `gdcm::Item::Write()`.


## 28.256 gdcmTransferSyntax.h File Reference

```
#include "gdcmSwapCode.h"
```

Include dependency graph for `gdcmTransferSyntax.h`:

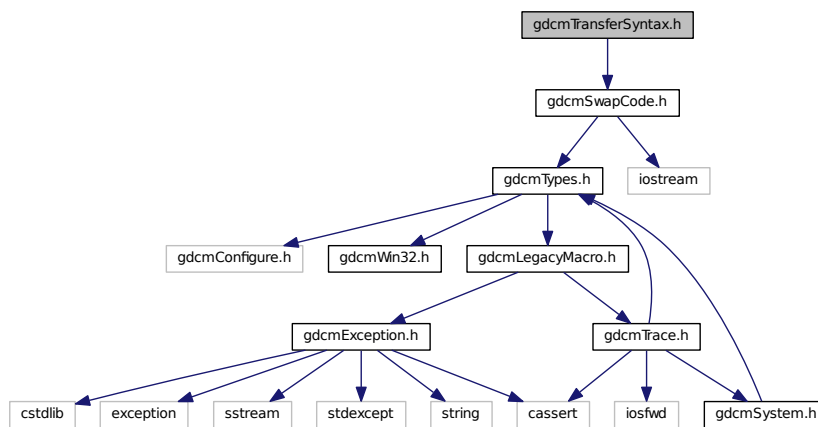

This graph shows which files directly or indirectly include this file:

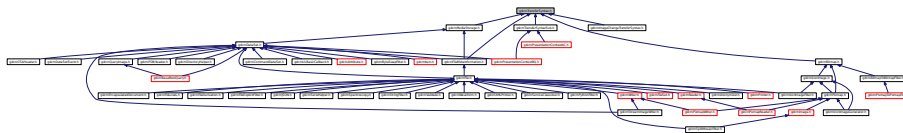

## Classes

- class [gdcm::TransferSyntax](#)

*Class to manipulate Transfer Syntax.*

## Namespaces


- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &_os, const TransferSyntax &ts)`

## 28.257 gdcmTransferSyntaxSub.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmTransferSyntax.h"
#include "gdcmUIDs.h"
Include dependency graph for gdcmTransferSyntaxSub.h:
```

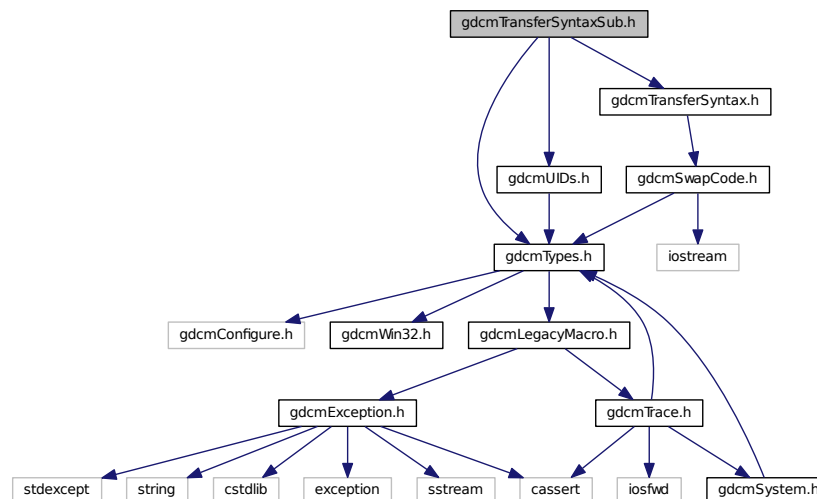

This graph shows which files directly or indirectly include this file:

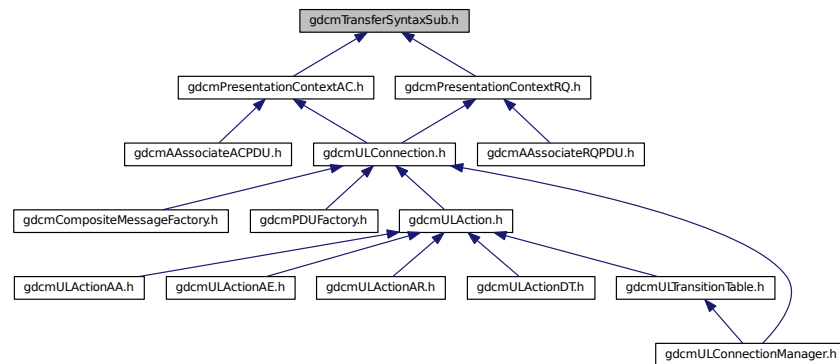

## Classes

- class [gdcm::network::TransferSyntaxSub](#)

*TransferSyntaxSub Table 9-15 TRANSFER SYNTAX SUB-ITEM FIELDS.*


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.258 gdcmType.h File Reference


```
#include "gdcmTypes.h"
#include <iostream>
```

Include dependency graph for gdcmType.h:

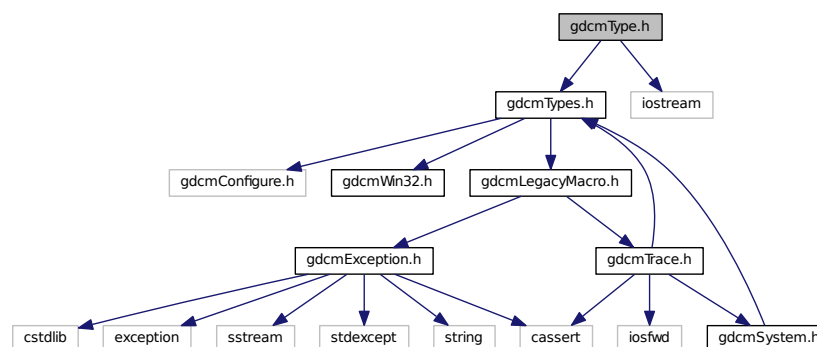




Include dependency graph for gdcmTypes.h:

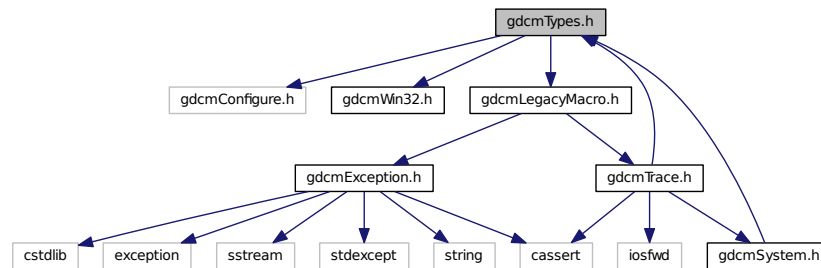

## 28.260 gdcmUIDGenerator.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmUIDGenerator.h:

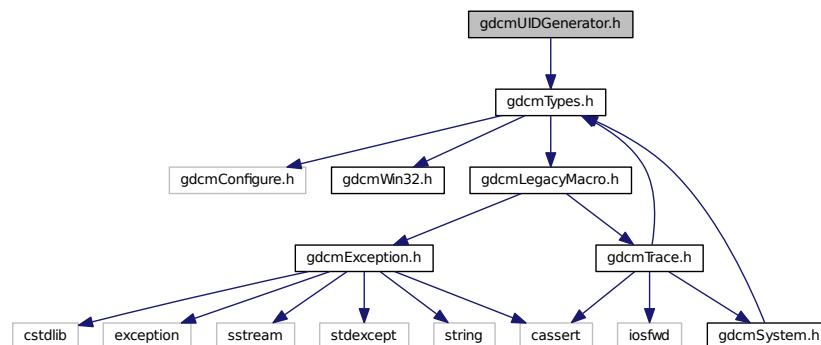

### Classes

- class [gdcm::UIDGenerator](#)  
*Class for generating unique UID.*


### Namespaces

- [gdcm](#)


## 28.261 gdcmUIDs.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for `gdcmUIDs.h`:

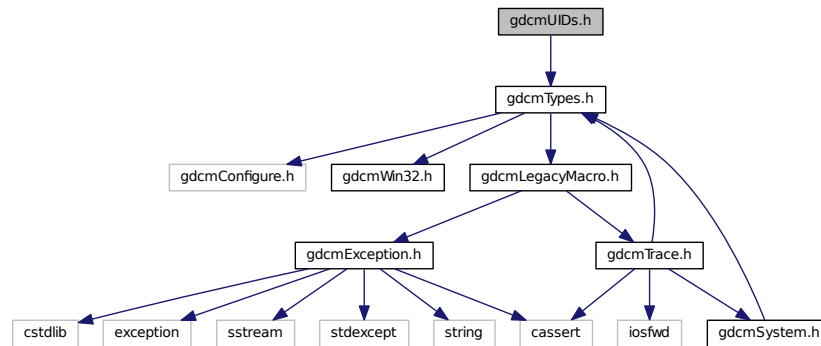

This graph shows which files directly or indirectly include this file:

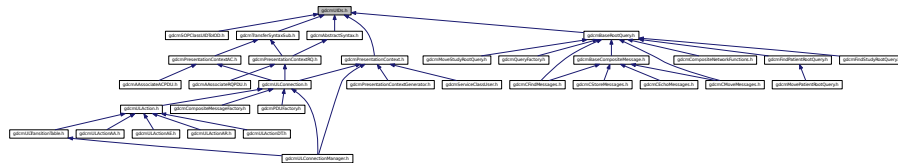

## Classes

- class `gdcm::UIDs`  
*all known uids*

## Namespaces

- `gdcm`

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const UIDs &uid)`


## 28.262 gdcmULAction.h File Reference

```


#include "gdcmNetworkStateID.h"
#include "gdcmULEvent.h"
#include "gdcmULConnection.h"

```

Include dependency graph for gdcmULAction.h:

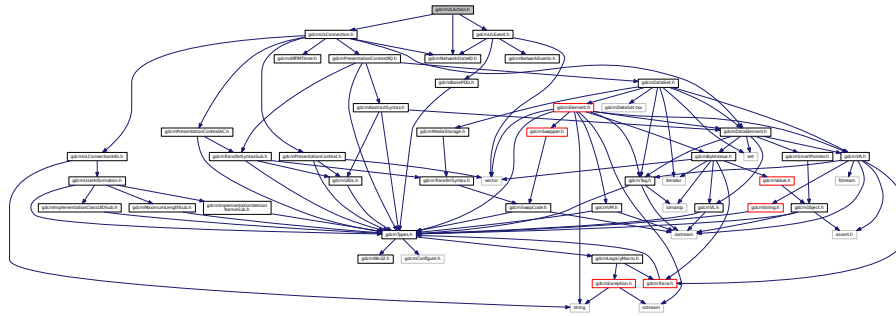

This graph shows which files directly or indirectly include this file:

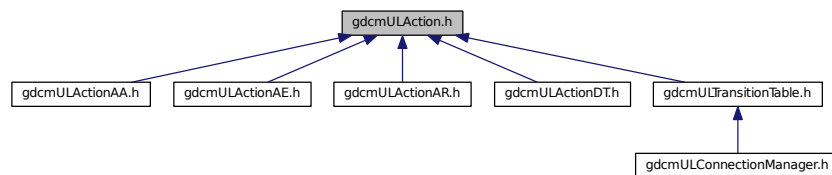

## Classes

- class [gdcm::network::ULAction](#)

*[ULAction](#) A [ULConnection](#) in a given [ULState](#) can perform certain [ULActions](#). This base class provides the interface for running those [ULActions](#) on a given [ULConnection](#).*


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.263 gdcmULActionAA.h File Reference

```
#include "gdcmULAction.h"
```

Include dependency graph for `gdcmlActionAA.h`:

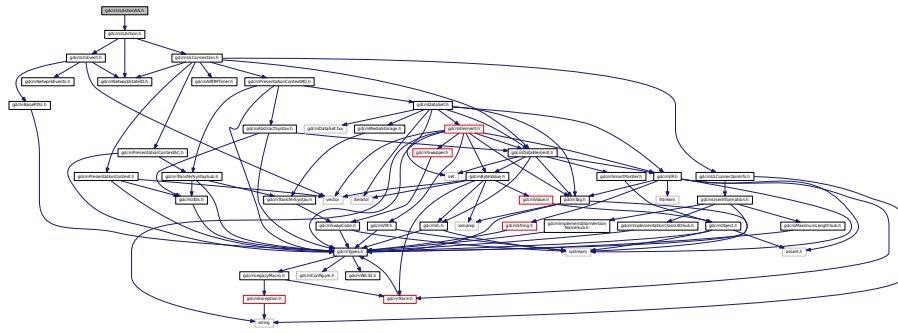

## Classes

- class `gdcml::network::ULActionAA1`
- class `gdcml::network::ULActionAA2`
- class `gdcml::network::ULActionAA3`
- class `gdcml::network::ULActionAA4`
- class `gdcml::network::ULActionAA5`
- class `gdcml::network::ULActionAA6`
- class `gdcml::network::ULActionAA7`
- class `gdcml::network::ULActionAA8`


## Namespaces

- `gdcml`
- `gdcml::network`

## 28.264 gdcmlActionAE.h File Reference

```
#include "gdcmlAction.h"
```

Include dependency graph for `gdcmlActionAE.h`:

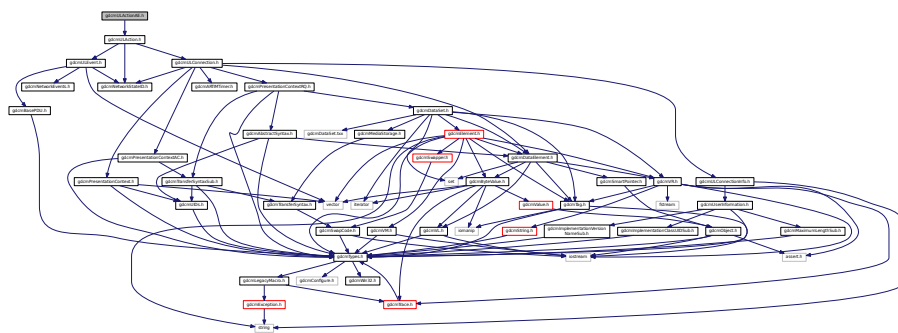



## Classes

- class [gdcm::network::ULActionAE1](#)
- class [gdcm::network::ULActionAE2](#)
- class [gdcm::network::ULActionAE3](#)
- class [gdcm::network::ULActionAE4](#)
- class [gdcm::network::ULActionAE5](#)
- class [gdcm::network::ULActionAE6](#)
- class [gdcm::network::ULActionAE7](#)
- class [gdcm::network::ULActionAE8](#)


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.265 gdcmULActionAR.h File Reference

```
#include "gdcmULAction.h"
```

Include dependency graph for gdcmULActionAR.h:

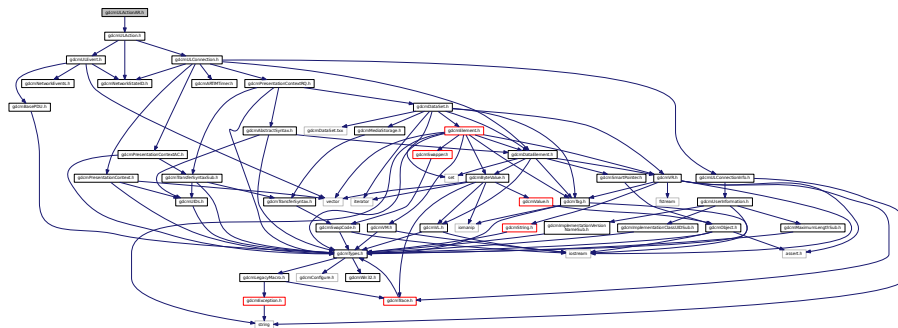

## Classes

- class [gdcm::network::ULActionAR1](#)
- class [gdcm::network::ULActionAR10](#)
- class [gdcm::network::ULActionAR2](#)
- class [gdcm::network::ULActionAR3](#)
- class [gdcm::network::ULActionAR4](#)
- class [gdcm::network::ULActionAR5](#)
- class [gdcm::network::ULActionAR6](#)
- class [gdcm::network::ULActionAR7](#)
- class [gdcm::network::ULActionAR8](#)
- class [gdcm::network::ULActionAR9](#)


## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.266 gdcmULActionDT.h File Reference

```
#include "gdcmULAction.h"
```

Include dependency graph for gdcmULActionDT.h:

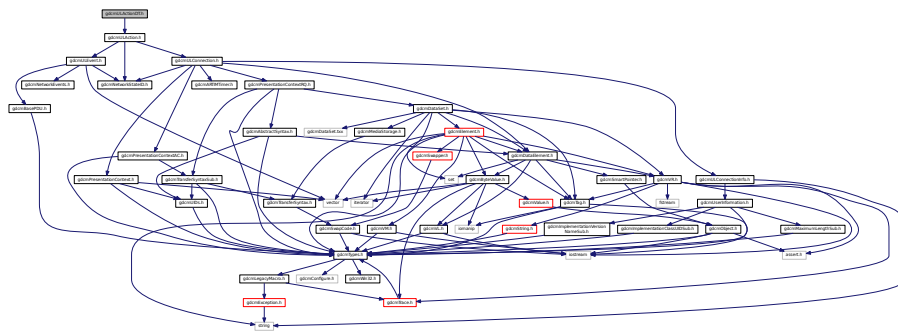

## Classes

- class [gdcm::network::ULActionDT1](#)
- class [gdcm::network::ULActionDT2](#)

## Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.267 gdcmULBasicCallback.h File Reference

```
#include "gdcmULConnectionCallback.h"
```

```
#include "gdcmDataSet.h"
```

```
#include <vector>
```

The diagram illustrates a complex dependency network for Glibc headers. At the top level, `glibcBasicCallback.h` serves as a primary header, branching into `glibcConnectorCallback.h` and `glibcDataSet.h`. These headers further develop into more specialized modules such as `glibcMediaStorage.h`, `glibcResource.h`, `glibcTransferSyntax.h`, `glibcSecurityCode.h`, `glibcSwapper.h`, `glibcTag.h`, `glibcValue.h`, `glibcTypes.h`, `glibcConfig.h`, `glibcMn32.h`, `glibcLgcnMacros.h`, `glibcException.h`, `glibcSystem.h`, and `glibcPace.h`. The network concludes with a series of lower-level or utility headers at the bottom, including `cstream`, `csong`, `stdexcept`, `cstdint`, `exception`, `casert`, `glibcSystem.h`, and `testad`. The dense web of arrows indicates extensive cross-dependencies throughout the entire system.


- class `gdcm::network::ULBasicCallback`

## Namespaces


- gdc
- gdc::network

```
#include "gdcNetworkStateID.h"
#include "gdcARTIMTimer.h"
#include "gdcULConnectionInfo.h"
#include "gdcPresentationContextRQ.h"
#include "gdcDataElement.h"
#include "gdcPresentationContextAC.h"
#include "gdcPresentationContext.h"
```

Include dependency graph for `gdcmULConnection.h`:

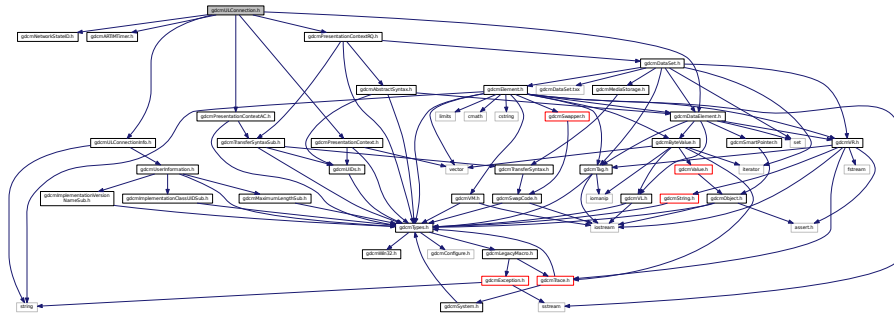

This graph shows which files directly or indirectly include this file:

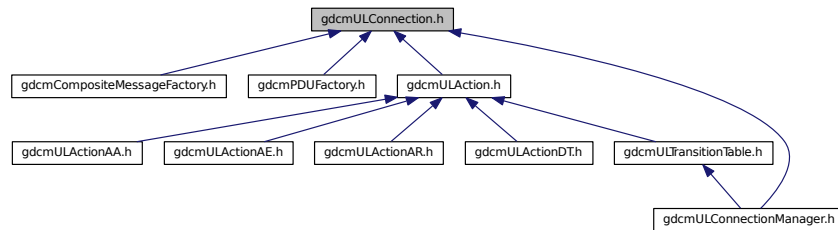

## Classes

- class `gdcm::network::ULConnection`

***ULConnection** This is the class that contains the socket to another machine, and passes data through itself, as well as maintaining a sense of state.*


## Namespaces

- `gdcm`
- `gdcm::network`


## 28.269 gdcmULConnectionCallback.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmULConnectionCallback.h:

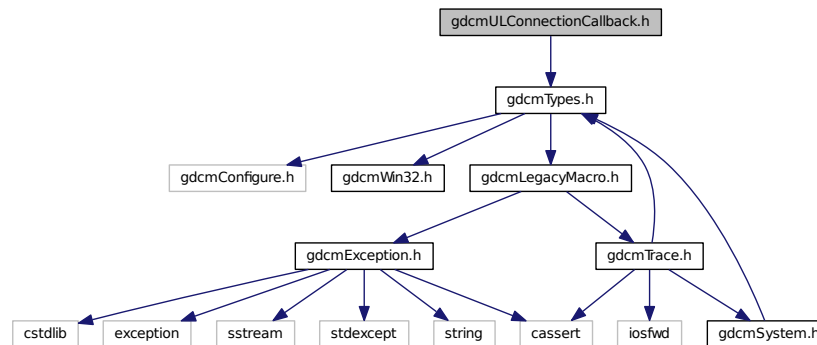

This graph shows which files directly or indirectly include this file:

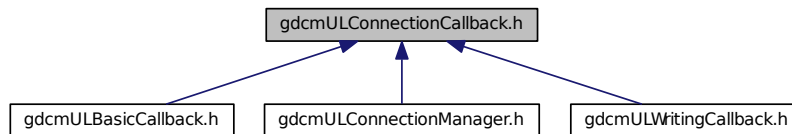

## Classes

- class [gdcm::network::ULConnectionCallback](#)


## Namespaces

- [gdcm](#)
- [gdcm::network](#)


## 28.270 gdcmULConnectionInfo.h File Reference

```
#include "gdcmUserInformation.h"
#include <string>
```

Include dependency graph for `gdcmULConnectionInfo.h`:

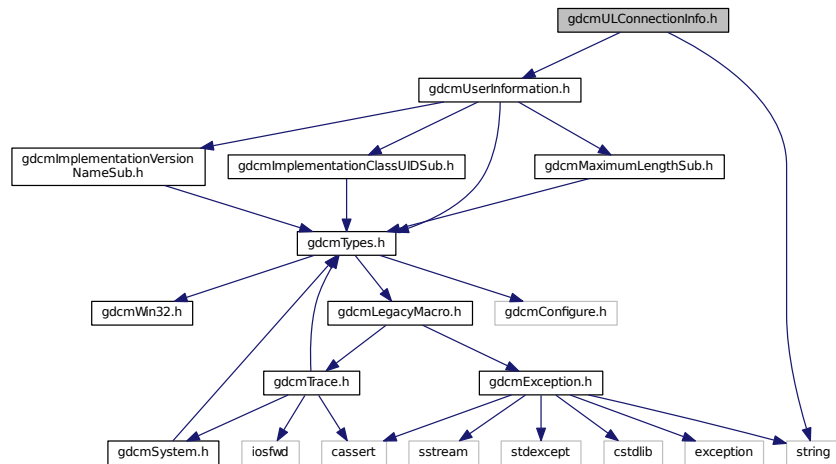

This graph shows which files directly or indirectly include this file:

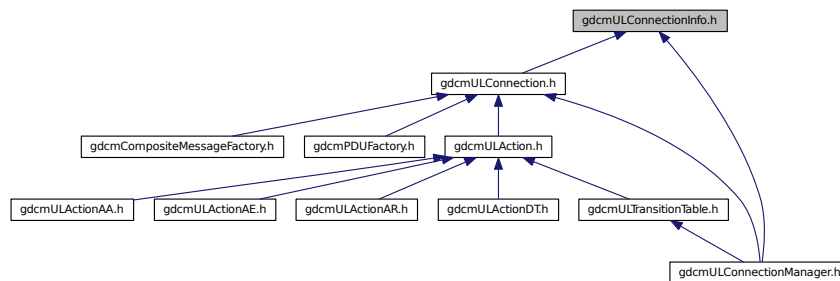

## Classes

- class `gdcm::network::ULConnectionInfo`


*`ULConnectionInfo` this class contains all the information about a particular connection as established by the user. That is, it's: User Information Calling AE Title Called AE Title IP address/computer name IP Port A connection must be established with this information, that's subsequently placed into various primitives for actual communication.*

## Namespaces

- `gdcm`
- `gdcm::network`

## 28.271 gdcmULConnectionManager.h File Reference

```
#include "gdcmULTransitionTable.h"
#include "gdcmULConnection.h"
#include "gdcmULConnectionInfo.h"
#include "gdcmPresentationDataValue.h"
#include "gdcmULConnectionCallback.h"
#include "gdcmSubject.h"
#include "gdcmPresentationContext.h"
Include dependency graph for gdcmULConnectionManager.h:
```

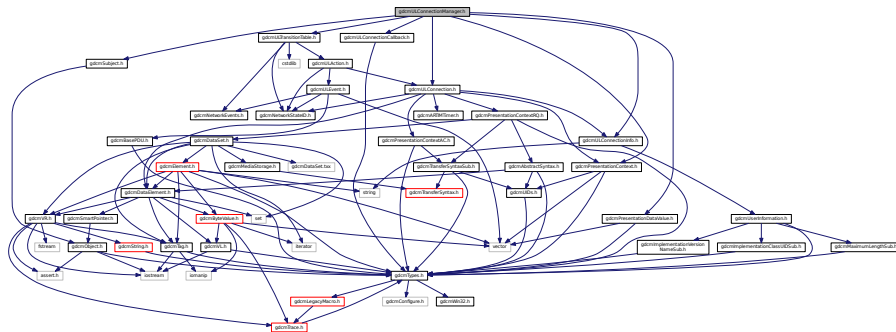

### Classes

- class [gdcm::network::ULConnectionManager](#)

*[ULConnectionManager](#) The [ULConnectionManager](#) performs actions on the [ULConnection](#) given inputs from the user and from the state of what's going on around the connection (ie, timeouts of the ARTIM timer, responses from the peer across the connection, etc).*


### Namespaces

- [gdcm](#)
- [gdcm::network](#)


## 28.272 gdcmULEvent.h File Reference

```
#include "gdcmNetworkStateID.h"
#include "gdcmNetworkEvents.h"
#include "gdcmBasePDU.h"
#include <vector>
```

Include dependency graph for `gdcmlEvent.h`:

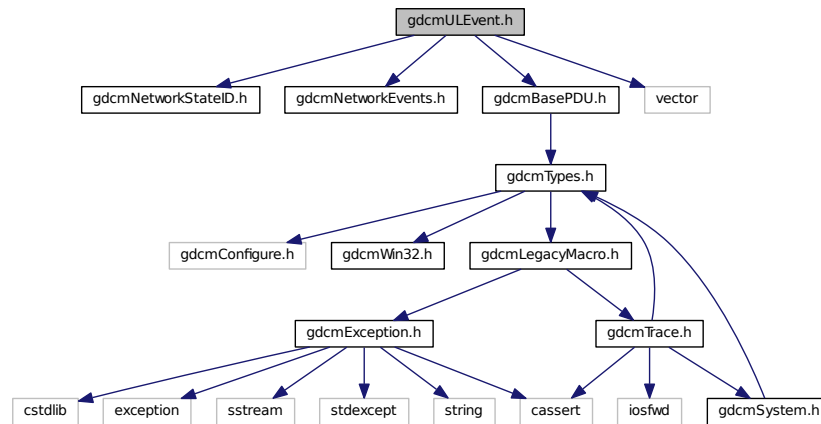

This graph shows which files directly or indirectly include this file:

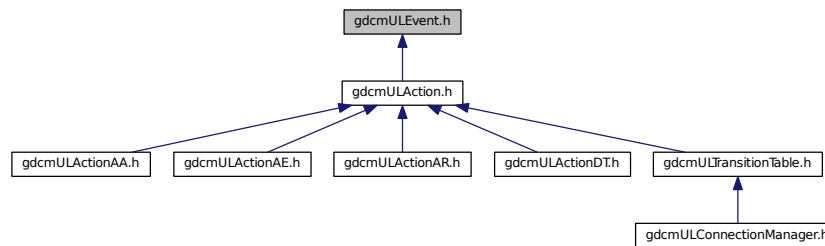

## Classes

- class `gdcml::network::ULEvent`  
*ULEvent* base class for network events.

## Namespaces

- `gdcml`
- `gdcml::network`


## 28.273 gdcmlULTransitionTable.h File Reference

```
#include "gdcmlNetworkStateID.h"
```




```
#include "gdcmNetworkEvents.h"
#include "gdcmULAction.h"
#include <cstdlib>
```

Include dependency graph for gdcmULTransitionTable.h:

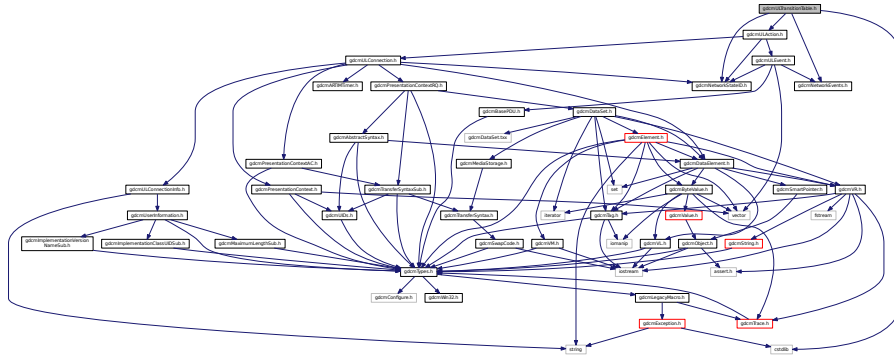

This graph shows which files directly or indirectly include this file:

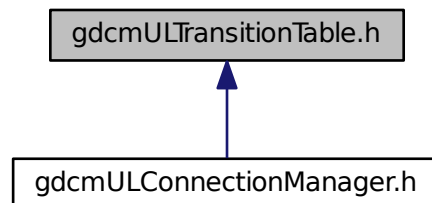

## Classes

- class [gdcm::network::TableRow](#)
- struct [gdcm::network::Transition](#)
- class [gdcm::network::ULTransitionTable](#)

*[ULTransitionTable](#) The transition table of all the ULEvents, new ULActions, and ULStates.*


## Namespaces

- [gdcm](#)
- [gdcm::network](#)





Include dependency graph for `gdcmUnpacker12Bits.h`:

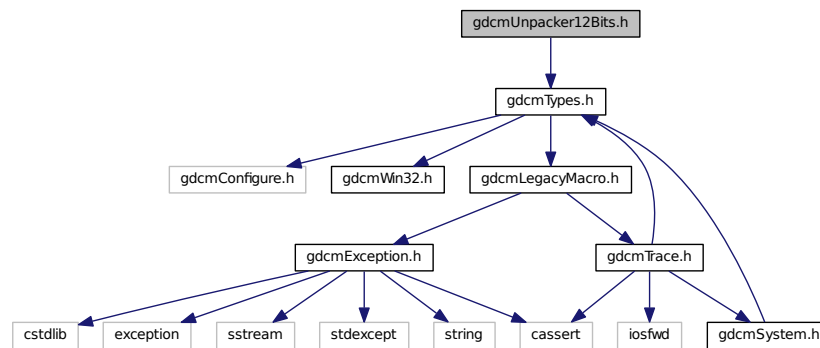

## Classes

- class `gdcm::Unpacker12Bits`  
*Pack/Unpack 12 bits pixel into 16bits.*

## Namespaces


- `gdcm`

## 28.278 gdcmUsage.h File Reference


```
#include "gdcmTypes.h"
```

```
#include <iostream>
```

Include dependency graph for `gdcmUsage.h`:

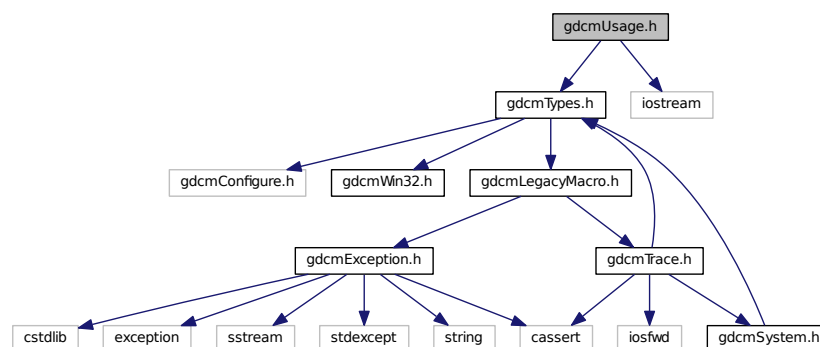

This graph shows which files directly or indirectly include this file:

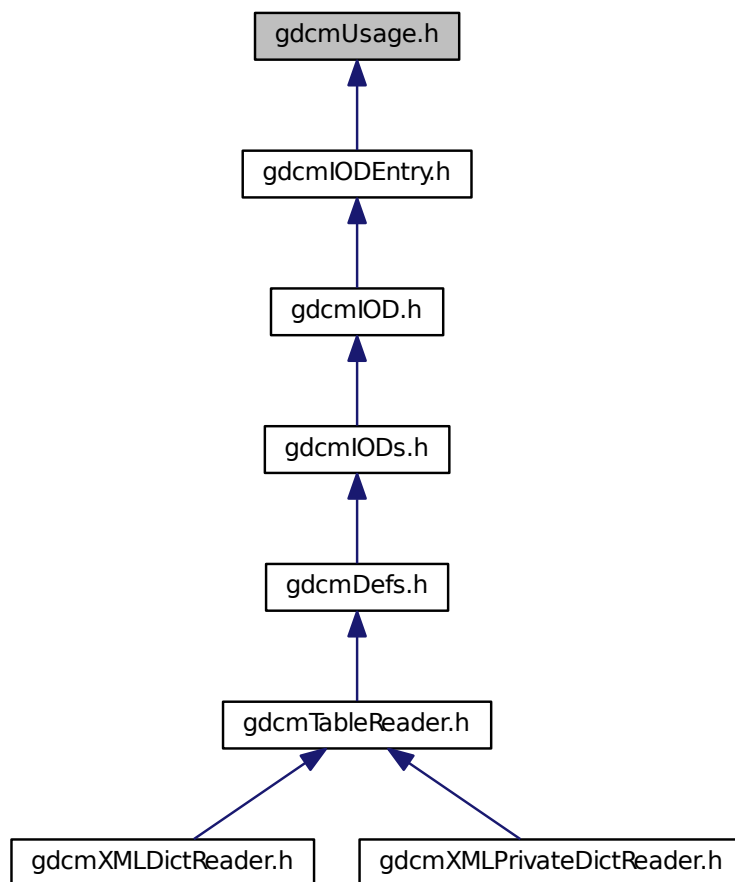

## Classes

- class `gdcm::Usage`  
*Usage.*

## Namespaces

- `gdcm`


## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const Usage &val)`


## 28.279 gdcmUserInformation.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmMaximumLengthSub.h"
#include "gdcmImplementationVersionNameSub.h"
#include "gdcmImplementationClassUIDSub.h"
```

Include dependency graph for gdcmUserInformation.h:

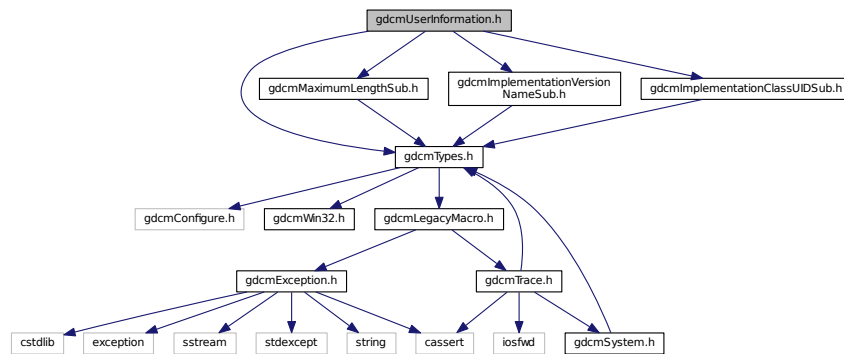

This graph shows which files directly or indirectly include this file:

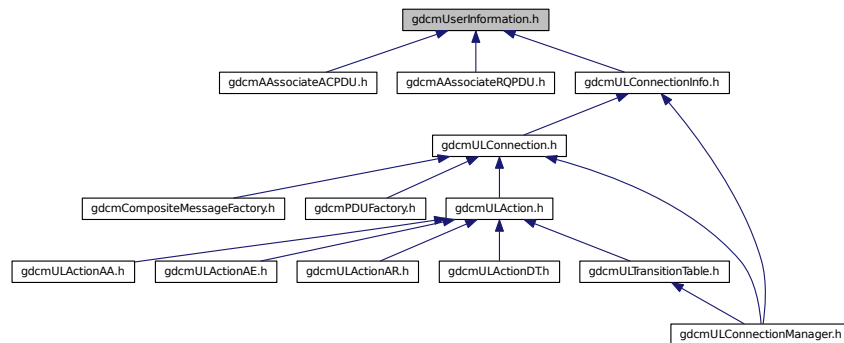

### Classes

- class [gdcm::network::UserInformation](#)  
*UserInformation Table 9-16 USER INFORMATION ITEM FIELDS.*


### Namespaces

- [gdcm](#)
- [gdcm::network](#)

## 28.280 gdcmUUIDGenerator.h File Reference

```
#include "gdcmTypes.h"
```

Include dependency graph for gdcmUUIDGenerator.h:

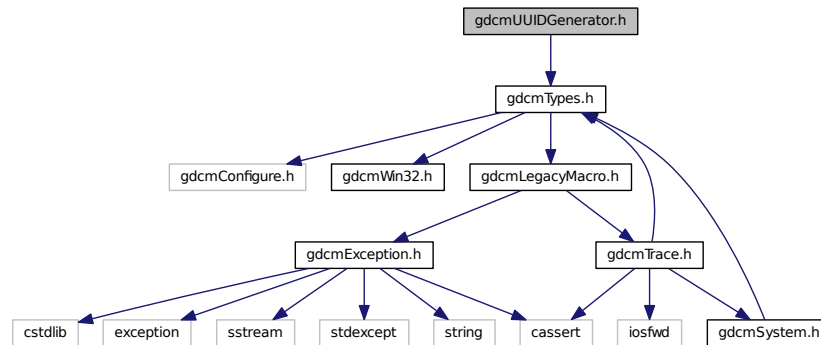

### Classes

- class `gdcm::UUIDGenerator`


*Class for generating unique UUID generate DCE 1.1 uid.*

### Namespaces

- `gdcm`

## 28.281 gdcmValidate.h File Reference


```
#include "gdcmFile.h"
```

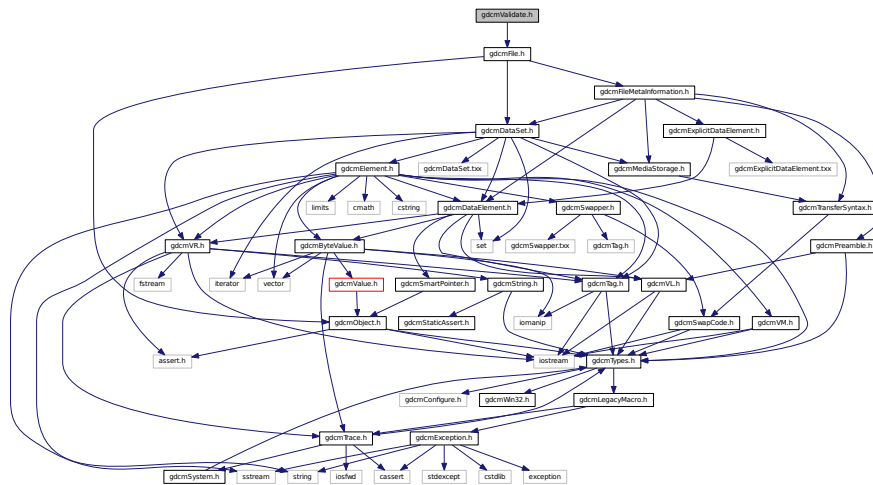


---



Include dependency graph for gdcmValue.h:

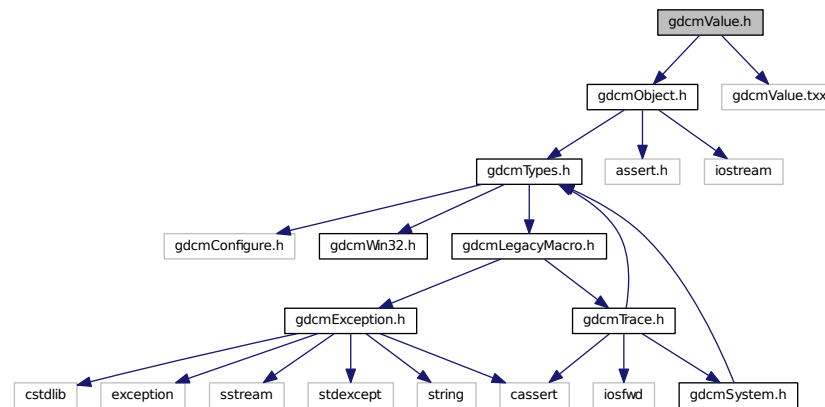

This graph shows which files directly or indirectly include this file:

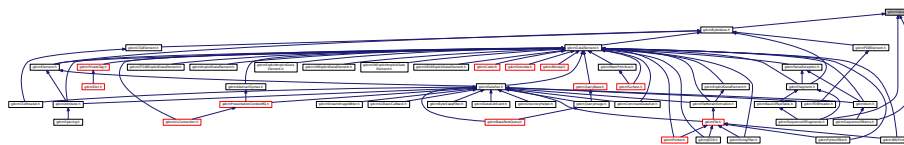

## Classes

- class [gdcm::Value](#)

*Class to represent the value of a Data [Element](#).*


## Namespaces

- [gdcm](#)

## 28.283 gdcmValueIO.h File Reference

```
#include "gdcmTypes.h"
#include "gdcmValueIO.txx"
```

Include dependency graph for `gdcmValueIO.h`:

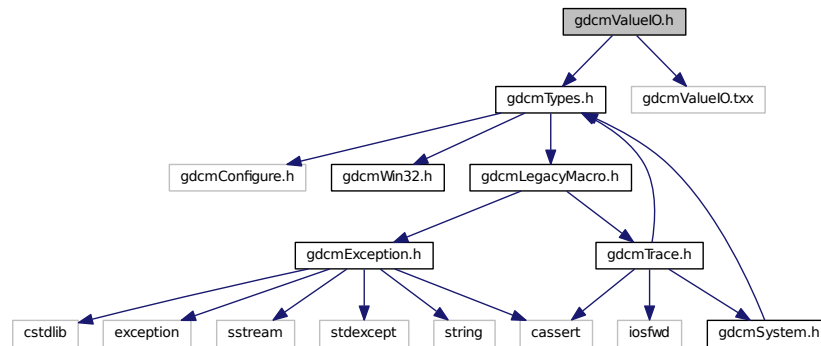

## Classes

- class [gdcm::ValueIO< TDE, TSwap, TType >](#)  
Class to dispatch template calls.

## Namespaces


- [gdcm](#)

## 28.284 gdcmVersion.h File Reference

```
#include "gdcmTypes.h"
```

```
#include <iostream>
```

Include dependency graph for `gdcmVersion.h`:

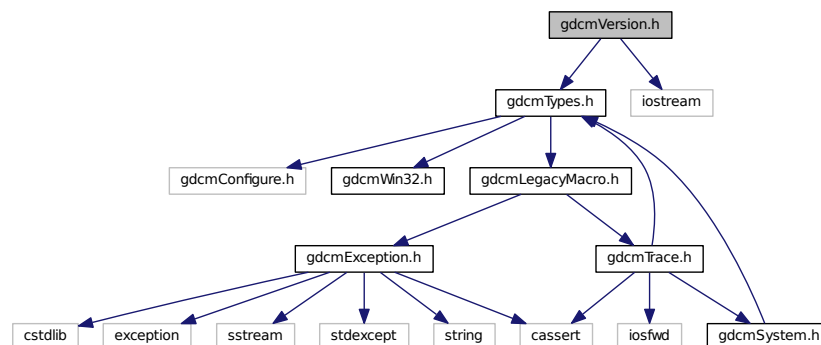

## Classes

- class [gdcm::Version](#)  
major/minor and build version

## Namespaces

- [gdcm](#)


## Functions

- std::ostream & [gdcm::operator<<](#) (std::ostream &os, const Version &v)


## 28.285 gdcviewer.dox File Reference

## 28.286 gdcmVL.h File Reference

```
#include "gdcmTypes.h"
#include <iostream>
Include dependency graph for gdcmVL.h:
```

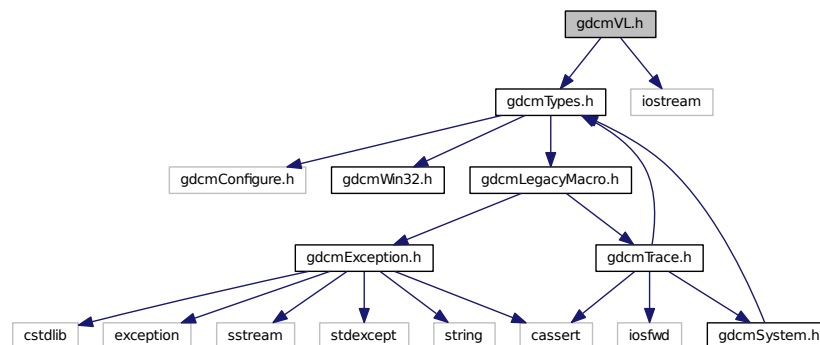

This graph shows which files directly or indirectly include this file:

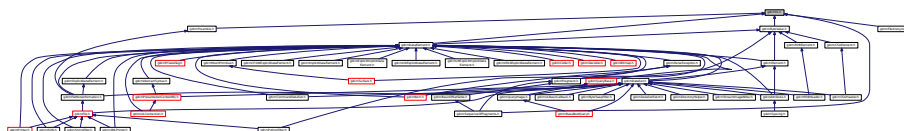

## Classes

- class [gdcm::VL](#)

*Value Length.*

## Namespaces

- [gdcm](#)

## Functions


- `std::ostream & gdcm::operator<< (std::ostream &os, const VL &val)`

## 28.287 gdcmVM.h File Reference


```
#include "gdcmTypes.h"
```

```
#include <iostream>
```

Include dependency graph for `gdcmVM.h`:

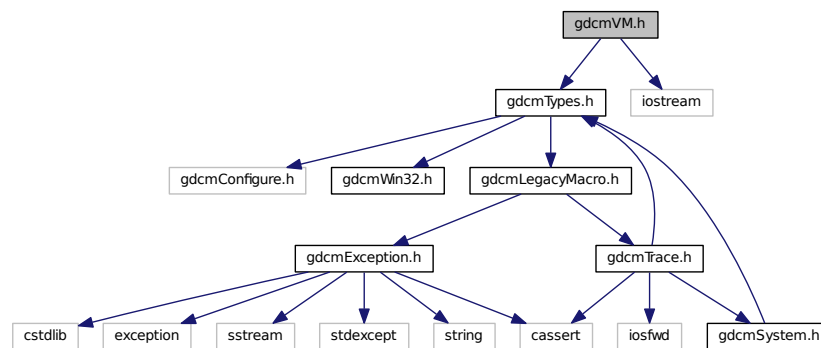

This graph shows which files directly or indirectly include this file:

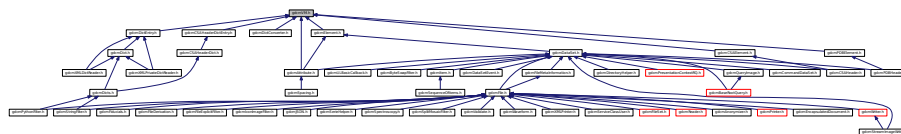

## Classes

- class [gdcm::VM](#)

*Value Multiplicity Looking at the DICOMV3 dict only there is very few cases: 1 2 3 4 5 6 8 16 24 1-2 1-3 1-8 1-32 1-99 1-n 2-2n 2-n 3-3n 3-n.*

- struct [gdcm::VMToLength< T >](#)

## Namespaces

- [gdcm](#)

## Macros

- `#define TYPETOLENGTH(type, length)`

## Functions

- `std::ostream & gdcm::operator<< (std::ostream &_os, const VM &_val)`

### 28.287.1 Macro Definition Documentation

#### 28.287.1.1 `#define TYPETOLENGTH( type, length )`


#### Value:

```
template<> struct VMToLength<VM::type> \
{ enum { Length = length }; };
```


## 28.288 gdcmVR.h File Reference

```
#include "gdcmTag.h"
#include "gdcmTrace.h"
#include "gdcmString.h"
#include <iostream>
#include <fstream>
#include <assert.h>
```

Include dependency graph for gdcmVR.h:

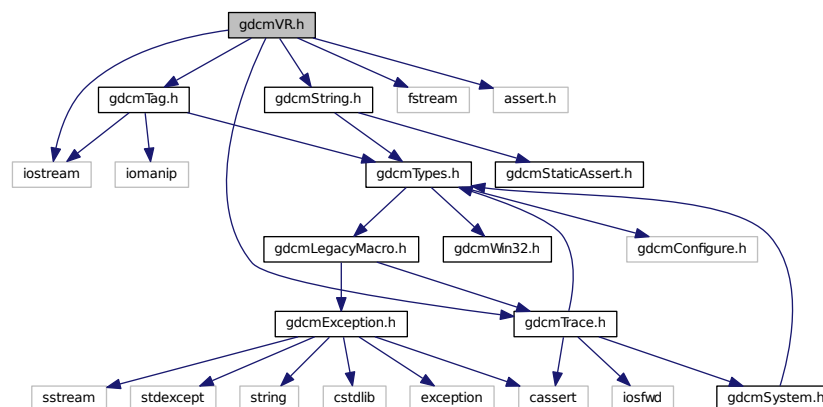

This graph shows which files directly or indirectly include this file:

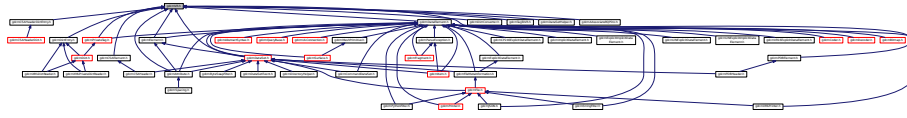

## Classes

- struct [gdcm::UI](#)
- class [gdcm::VR](#)

*VR class This is adapted from DICOM standard The biggest difference is the INVALID VR and the composite one that differ from standard (more like an addition) This allow us to represent all the possible case express in the DICOMV3 dict.*

- struct [gdcm::VRToEncoding< T >](#)
- struct [gdcm::VRToType< T >](#)

## Namespaces

- [gdcm](#)

## Macros

- #define [TYPETOENCODING](#)(type, rep, rtype)
- #define [VRTypeTemplateCase](#)(type)

## Typedefs

- typedef String<"\\", 16 > [gdcm::AECComp](#)
- typedef String<"\\", 64 > [gdcm::ASComp](#)
- typedef String<"\\", 16 > [gdcm::CSCComp](#)
- typedef String<"\\", 64 > [gdcm::DACComp](#)
- typedef String<"\\", 64 > [gdcm::DTComp](#)
- typedef String<"\\", 64 > [gdcm::LOComp](#)
- typedef String<"\\", 64 > [gdcm::LTComp](#)
- typedef String<"\\", 64 > [gdcm::PNComp](#)
- typedef String<"\\", 64 > [gdcm::SHComp](#)
- typedef String<"\\", 64 > [gdcm::STComp](#)
- typedef String<"\\", 16 > [gdcm::TMComp](#)
- typedef String<"\\", 64, 0 > [gdcm::UIComp](#)
- typedef String<"\\", 64 > [gdcm::UTComp](#)

## Functions

- std::ostream & [gdcm::operator<<](#) (std::ostream &\_os, const VR &val)
- std::ostream & [gdcm::operator<<](#) (std::ostream &\_os, const UI &\_val)
- [gdcm::TYPETOENCODING](#) (SQ, VRBINARY, unsigned char) TYPETOENCODING(UN









## 28.292 gdcmWriter.h File Reference

```
#include "gdcmFile.h"
```

Include dependency graph for gdcmWriter.h:

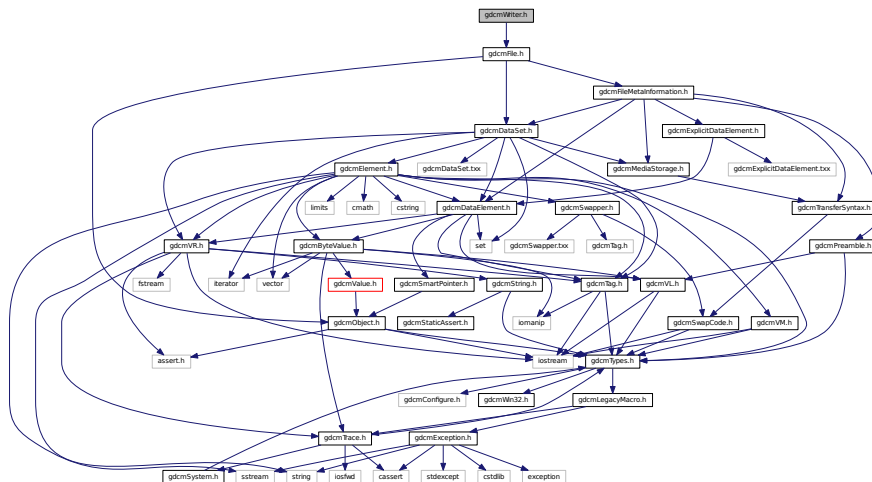

This graph shows which files directly or indirectly include this file:

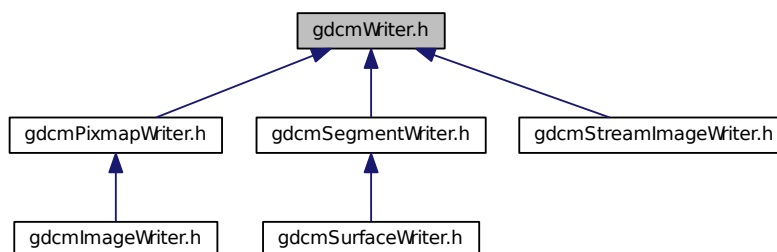

## Classes

- class [gdcm::Writer](#)


*Writer* ala DOM (Document *Object* Model) This class is a non-validating writer, it will only performs well- formedness check only.

## Namespaces

- [gdcm](#)



Include dependency graph for gdcmXMLPrinter.h:

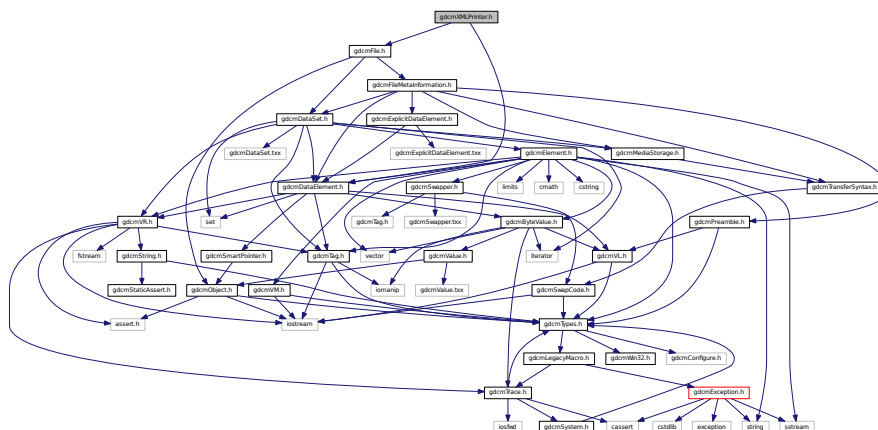

## Classes

- class [gdcm::XMLPrinter](#)


## Namespaces

- [gdcm](#)

## 28.296 gdcmXMLPrivateDictReader.h File Reference

```
#include "gdcmTableReader.h"
#include "gdcmDict.h"
#include "gdcmDictEntry.h"
#include "gdcmTag.h"
```

Include dependency graph for gdcmXMLPrivateDictReader.h:

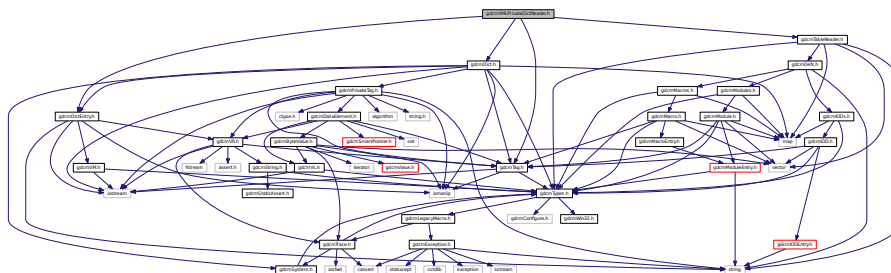

## Classes

- class [gdcm::XMLPrivateDictReader](#)  
Class for representing a *XMLPrivateDictReader*.

## Namespaces


- [gdcm](#)

## 28.297 README.txt File Reference


## 28.298 TestsList.txt File Reference

## 28.299 vtkGDCMImageReader.h File Reference

```
#include "vtkMedicalImageReader2.h"  
#include "vtkImageData.h"  
Include dependency graph for vtkGDCMImageReader.h:
```

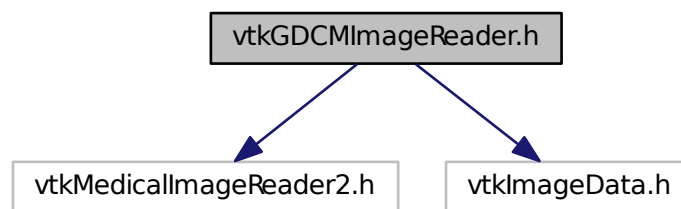

This graph shows which files directly or indirectly include this file:

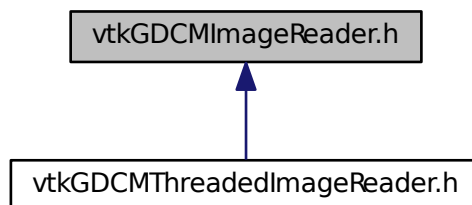

## Classes

- class [vtkGDCMImageReader](#)

## Namespaces

- [gdcm](#)

## Macros

- `#define VTK_CMYK` 8
- `#define VTK_INVERSE_LUMINANCE` 5
- `#define VTK_LOOKUP_TABLE` 6
- `#define VTK_YBR` 7

### 28.299.1 Macro Definition Documentation

28.299.1.1 `#define VTK_CMYK` 8

28.299.1.2 `#define VTK_INVERSE_LUMINANCE` 5

28.299.1.3 `#define VTK_LOOKUP_TABLE` 6


28.299.1.4 `#define VTK_YBR` 7

## 28.300 vtkGDCMImageReader2.h File Reference

```
#include "vtkMedicalImageReader2.h"
```

```
#include "vtkImageData.h"
```

Include dependency graph for vtkGDCMImageReader2.h:

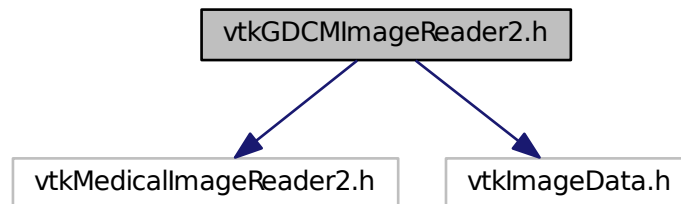

## Classes

- class [vtkGDCMImageReader2](#)

## Namespaces

- [gdcm](#)

## Macros

- `#define VTK_CMYK 8`
- `#define VTK_INVERSE_LUMINANCE 5`
- `#define VTK_LOOKUP_TABLE 6`
- `#define VTK_YBR 7`

### 28.300.1 Macro Definition Documentation

28.300.1.1 `#define VTK_CMYK 8`

28.300.1.2 `#define VTK_INVERSE_LUMINANCE 5`


28.300.1.3 `#define VTK_LOOKUP_TABLE 6`

28.300.1.4 `#define VTK_YBR 7`

## 28.301 vtkGDCMImageWriter.h File Reference

```
#include "vtkImageWriter.h"
```

Include dependency graph for vtkGDCMImageWriter.h:

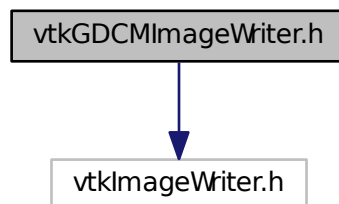


## Classes

- class `vtkGDCMImageWriter`

## 28.302 vtkGDCMMedicalImageProperties.h File Reference

```
#include "vtkMedicalImageProperties.h"
```

Include dependency graph for vtkGDCMMedicalImageProperties.h:

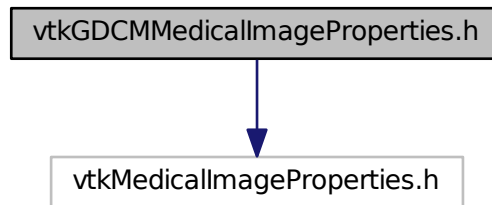

### Classes

- class [vtkGDCMMedicalImageProperties](#)


### Namespaces

- [gdc](#)

## 28.303 vtkGDCMPolyDataReader.h File Reference

```
#include "vtkPolyDataAlgorithm.h"
```

Include dependency graph for vtkGDCMPolyDataReader.h:

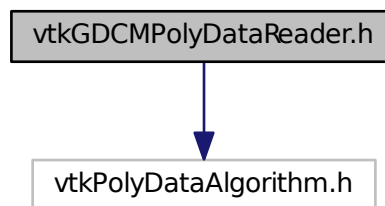

### Classes


- class [vtkGDCMPolyDataReader](#)

## Namespaces

- [gdcm](#)

## 28.304 vtkGDCMPolyDataWriter.h File Reference

```
#include "vtkPolyDataWriter.h"  
#include "vtkStringArray.h"  
#include "vtkStdString.h"  
Include dependency graph for vtkGDCMPolyDataWriter.h:
```

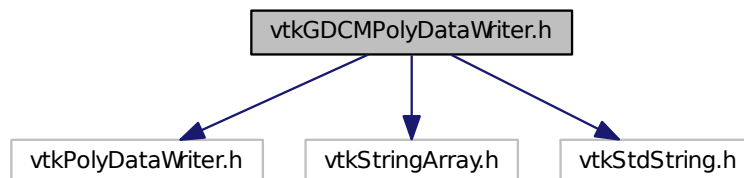

## Classes

- class [vtkGDCMPolyDataWriter](#)

## Namespaces


- [gdcm](#)

## 28.305 vtkGDCMTesting.h File Reference

```
#include "vtkObject.h"
```



Include dependency graph for vtkGDCMTesting.h:

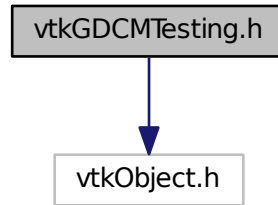


### Classes

- class [vtkGDCMTesting](#)

## 28.306 vtkGDCMThreadedImageReader.h File Reference

```
#include "vtkGDCMImageReader.h"
```

Include dependency graph for vtkGDCMThreadedImageReader.h:

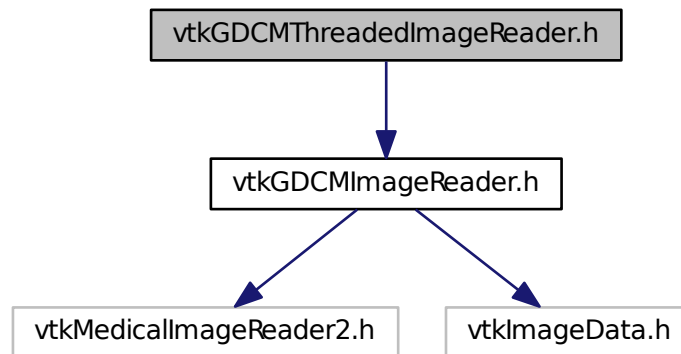


### Classes

- class [vtkGDCMThreadedImageReader](#)

## 28.307 vtkGDCMThreadedImageReader2.h File Reference

```
#include "vtkThreadedImageAlgorithm.h"
```

Include dependency graph for vtkGDCMThreadedImageReader2.h:

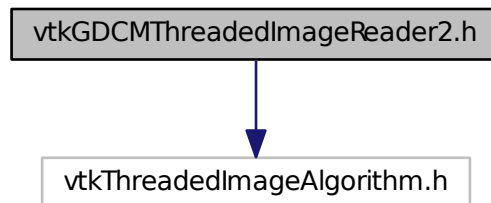


### Classes

- class [vtkGDCMThreadedImageReader2](#)

## 28.308 vtkImageColorViewer.h File Reference

```
#include "vtkObject.h"
```

Include dependency graph for vtkImageColorViewer.h:

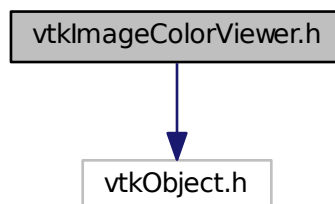


### Classes

- class [vtkImageColorViewer](#)

## 28.309 vtkImageMapToColors16.h File Reference

```
#include "vtkThreadedImageAlgorithm.h"
```

Include dependency graph for vtkImageMapToColors16.h:

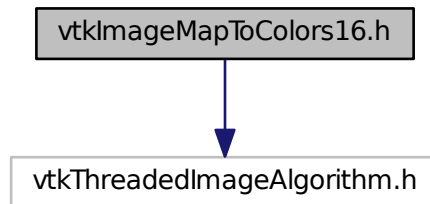


### Classes

- class [vtkImageMapToColors16](#)

## 28.310 vtkImageMapToWindowLevelColors2.h File Reference

```
#include "vtkImageMapToColors.h"
```

Include dependency graph for vtkImageMapToWindowLevelColors2.h:

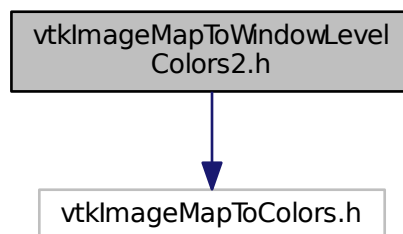


### Classes

- class [vtkImageMapToWindowLevelColors2](#)

## 28.311 vtkImagePlanarComponentsToComponents.h File Reference

```
#include "vtkImageAlgorithm.h"
```

Include dependency graph for vtkImagePlanarComponentsToComponents.h:

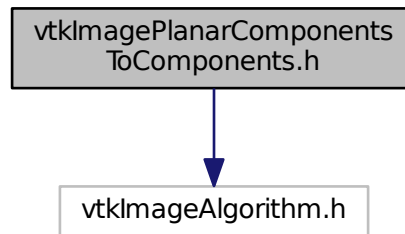


### Classes

- class [vtkImagePlanarComponentsToComponents](#)

## 28.312 vtkImageRGBToYBR.h File Reference

```
#include "vtkThreadedImageAlgorithm.h"
```

Include dependency graph for vtkImageRGBToYBR.h:

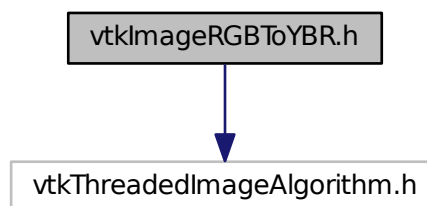


### Classes

- class [vtkImageRGBToYBR](#)

## 28.313 vtkImageYBRTToRGB.h File Reference

```
#include "vtkThreadedImageAlgorithm.h"
```

Include dependency graph for vtkImageYBRTToRGB.h:

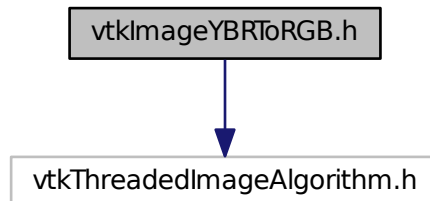

### Classes


- class [vtkImageYBRTToRGB](#)

## 28.314 vtkLookupTable16.h File Reference

```
#include "vtkLookupTable.h"
```

```
#include "vtkUnsignedShortArray.h"
```

Include dependency graph for vtkLookupTable16.h:

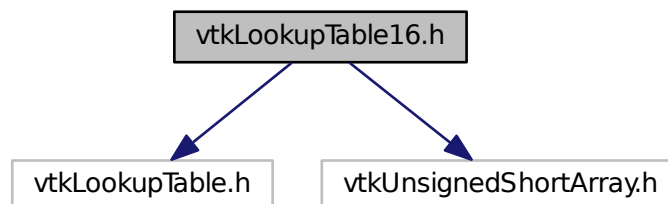


### Classes

- class [vtkLookupTable16](#)

## 28.315 vtkRTStructSetProperties.h File Reference

```
#include "vtkObject.h"
```

Include dependency graph for vtkRTStructSetProperties.h:

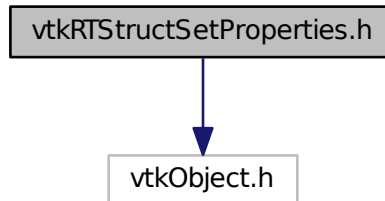

### Classes

- class [vtkRTStructSetProperties](#)

## Chapter 29

# Example Documentation

### 29.1 AWTMedical3.java

```
/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
package examples;

import vtk.*;
//import gdcm.*;

import vtk.util.VtkPanelContainer;
import vtk.util.VtkPanelUtil;
import vtk.util.VtkUtil;

import java.util.ArrayList;

import javax.swing.*;
import java.awt.*;
import java.io.File;

public class AWTMedical3 extends JComponent implements VtkPanelContainer {

    private vtkPanel renWin;

    vtkImageData ReadDataFile(File inSelectedFile){

        vtkImageData outImageData = null;
        Directory theDir = new Directory();

        String theInputDirectory = inSelectedFile.getPath();
        theDir.Load(theInputDirectory);

        Scanner theScanner = new Scanner();
        Tag theStudyTag = new Tag(0x0020,0x000d);
        Tag theSeriesTag = new Tag(0x0020,0x000e);
        theScanner.AddTag(theStudyTag); //get studies,
        theScanner.AddTag(theSeriesTag); //get studies,
        theScanner.Scan(theDir.GetFilenames());

        FilenamesType theStudyValues = theScanner.GetOrderedValues(theStudyTag);
        long theNumStudies = theStudyValues.size();
        //for now, take the first study, and nothing else.
        //and the return is actually not FilenamesType, just a
        //vector of strings
    }
}
```

```

    if (theNumStudies != 1)
        return outImageData;
    String theStudyVal = theStudyValues.get(0);
    //now, get all the values from the scanner that are in that
    //study, then from that get their different series
    FilenamesType theFilenames =
        theScanner.GetAllFilenamesFromTagToValue(theStudyTag, theStudyVal);

    //from that set of filenames, isolate individual series
    //conclude that singleton series = RT struct (can do further
    //checking for things like MIPs and the like)
    //and multiple series entries = volumetric data
    theScanner.Scan(theFilenames);
    FilenamesType theSeriesValues = theScanner.GetOrderedValues(theSeriesTag);
    String studyUID = theScanner.GetValue(theScanner.GetFilenames().get(0), theStudyTag);
    long theNumSeries = theSeriesValues.size();
    for (int i = 0; i < theNumSeries; i++) {
        FilenamesType theSeriesFiles =
            theScanner.GetAllFilenamesFromTagToValue(theSeriesTag, theSeriesValues.get(i));
        long theNumFilesInSeries = theSeriesFiles.size();
        if (theNumFilesInSeries > 1) { //assume it's CT or volumetric data
            //for now, assume a single volume
            //could have multiples, like PET and CT

            IPPSorter sorter = new IPPSorter();
            sorter.SetComputeZSpacing(true);
            sorter.SetZSpacingTolerance(0.001);
            Boolean sorted = sorter.Sort(theSeriesFiles);
            if (!sorted){
                //need some better way to handle failures here
                return outImageData;
            }

            FilenamesType sortedFT = sorter.GetFilenames();
            long theSize = sortedFT.size();
            vtkStringArray sa = new vtkStringArray();
            ArrayList<String> theStrings = new ArrayList<String>();

            vtkGDCMImageReader gdcmReader = new
            vtkGDCMImageReader();
            for (int j = 0; j < theSize; j++) {
                String theFileName = sortedFT.get(j);
                if (gdcmReader.CanReadFile(theFileName) > 0){
                    theStrings.add(theFileName);
                    sa.InsertNextValue(theFileName);
                } else {
                    //this is a busted series
                    //need some more appropriate error here
                    return outImageData;
                }
            }

            gdcmReader.SetFileNames(sa);

            gdcmReader.Update();

            outImageData = gdcmReader.GetOutput(); //the zeroth output should be the image
        }
    }
    String theImageInfo = "";
    if (outImageData != null){
        theImageInfo = outImageData.Print();
    }
    return outImageData;
}

//this function is a rewrite of Medical3 to see if data can
//be loaded via gdcm easily
public AWTMedical3(File inFile) {
    // Create the buttons.
    renWin = new vtkPanel();

    vtkImageData theImageData = ReadDataFile(inFile);

    // An isosurface, or contour value of 500 is known to correspond to the
    // skin of the patient. Once generated, a vtkPolyDataNormals filter is
    // is used to create normals for smooth surface shading during rendering.
    // The triangle stripper is used to create triangle strips from the
    // isosurface these render much faster on some systems.
    vtkContourFilter skinExtractor = new vtkContourFilter();
    skinExtractor.SetInput(theImageData);

```



```

skinExtractor.SetValue(0, 500);
vtkPolyDataNormals skinNormals = new vtkPolyDataNormals();
skinNormals.SetInput(skinExtractor.GetOutput());
skinNormals.SetFeatureAngle(60.0);
//      vtkStripper skinStripper = new vtkStripper();
//      skinStripper.SetInput(skinNormals.GetOutput());
vtkPolyDataMapper skinMapper = new vtkPolyDataMapper();
skinMapper.SetInput(skinNormals.GetOutput());
skinMapper.ScalarVisibilityOff();
vtkActor skin = new vtkActor();
skin.SetMapper(skinMapper);
skin.GetProperty().SetDiffuseColor(1, .49, .25);
skin.GetProperty().SetSpecular(.3);
skin.GetProperty().SetSpecularPower(20);

// An isosurface, or contour value of 1150 is known to correspond to the
// skin of the patient. Once generated, a vtkPolyDataNormals filter is
// is used to create normals for smooth surface shading during rendering.
// The triangle stripper is used to create triangle strips from the
// isosurface these render much faster on some systems.
vtkContourFilter boneExtractor = new vtkContourFilter();
boneExtractor.SetInput(theImageData);
boneExtractor.SetValue(0, 1150);
vtkPolyDataNormals boneNormals = new vtkPolyDataNormals();
boneNormals.SetInput(boneExtractor.GetOutput());
boneNormals.SetFeatureAngle(60.0);
vtkStripper boneStripper = new vtkStripper();
boneStripper.SetInput(boneNormals.GetOutput());
vtkPolyDataMapper boneMapper = new vtkPolyDataMapper();
boneMapper.SetInput(boneStripper.GetOutput());
boneMapper.ScalarVisibilityOff();
vtkActor bone = new vtkActor();
bone.SetMapper(boneMapper);
bone.GetProperty().SetDiffuseColor(1, 1, .9412);

// An outline provides context around the data.
vtkOutlineFilter outlineData = new vtkOutlineFilter();
outlineData.SetInput(theImageData);
vtkPolyDataMapper mapOutline = new vtkPolyDataMapper();
mapOutline.SetInput(outlineData.GetOutput());
vtkActor outline = new vtkActor();
outline.SetMapper(mapOutline);
outline.GetProperty().SetColor(0, 0, 0);

// Now we are creating three orthogonal planes passing through the
// volume. Each plane uses a different texture map and therefore has
// different coloration.

// Start by creating a black/white lookup table.
vtkLookupTable bwLut = new vtkLookupTable();
bwLut.SetTableRange(0, 2000);
bwLut.SetSaturationRange(0, 0);
bwLut.SetHueRange(0, 0);
bwLut.SetValueRange(0, 1);
bwLut.Build();

// Now create a lookup table that consists of the full hue circle (from
// HSV);.
vtkLookupTable hueLut = new vtkLookupTable();
hueLut.SetTableRange(0, 2000);
hueLut.SetHueRange(0, 1);
hueLut.SetSaturationRange(1, 1);
hueLut.SetValueRange(1, 1);
hueLut.Build();

// Finally, create a lookup table with a single hue but having a range
// in the saturation of the hue.
vtkLookupTable satLut = new vtkLookupTable();
satLut.SetTableRange(0, 2000);
satLut.SetHueRange(.6, .6);
satLut.SetSaturationRange(0, 1);
satLut.SetValueRange(1, 1);
satLut.Build();

// Create the first of the three planes. The filter vtkImageMapToColors
// maps the data through the corresponding lookup table created above.
// The vtkImageActor is a type of vtkProp and conveniently displays an
// image on a single quadrilateral plane. It does this using texture
// mapping and as a result is quite fast. (Note: the input image has to
// be unsigned char values, which the vtkImageMapToColors produces.);
// Note also that by specifying the DisplayExtent, the pipeline

```

```

// requests data of this extent and the vtkImageMapToColors only
// processes a slice of data.
vtkImageMapToColors sagittalColors = new vtkImageMapToColors();
sagittalColors.SetInput(theImageData);
sagittalColors.SetLookupTable(bwLut);
vtkImageActor sagittal = new vtkImageActor();
sagittal.SetInput(sagittalColors.GetOutput());
sagittal.SetDisplayExtent(32, 32, 0, 63, 0, 92);

// Create the second (axial); plane of the three planes. We use the same
// approach as before except that the extent differs.
vtkImageMapToColors axialColors = new vtkImageMapToColors();
axialColors.SetInput(theImageData);
axialColors.SetLookupTable(hueLut);
vtkImageActor axial = new vtkImageActor();
axial.SetInput(axialColors.GetOutput());
axial.SetDisplayExtent(0, 63, 0, 63, 46, 46);

// Create the third (coronal); plane of the three planes. We use the same
// approach as before except that the extent differs.
vtkImageMapToColors coronalColors = new vtkImageMapToColors();
coronalColors.SetInput(theImageData);
coronalColors.SetLookupTable(satLut);
vtkImageActor coronal = new vtkImageActor();
coronal.SetInput(coronalColors.GetOutput());
coronal.SetDisplayExtent(0, 63, 32, 32, 0, 92);

// It is convenient to create an initial view of the data. The FocalPoint
// and Position form a vector direction. Later on (ResetCamera() method)
// this vector is used to position the camera to look at the data in
// this direction.
vtkCamera aCamera = new vtkCamera();
aCamera.SetViewUp(0, 0, -1);
aCamera.SetPosition(0, 1, 0);
aCamera.SetFocalPoint(0, 0, 0);
aCamera.ComputeViewPlaneNormal();

// Actors are added to the renderer. An initial camera view is created.
// The Dolly() method moves the camera towards the FocalPoint,
// thereby enlarging the image.
renWin.GetRenderer().AddActor(sagittal);
renWin.GetRenderer().AddActor(axial);
renWin.GetRenderer().AddActor(coronal);
renWin.GetRenderer().AddActor(outline);
renWin.GetRenderer().AddActor(skin);
renWin.GetRenderer().AddActor(bone);

// Turn off bone for this example.
bone.VisibilityOff();

// Set skin to semi-transparent.
skin.GetProperty().SetOpacity(0.5);

// An initial camera view is created. The Dolly() method moves
// the camera towards the FocalPoint, thereby enlarging the image.
renWin.GetRenderer().SetActiveCamera(aCamera);
renWin.GetRenderer().ResetCamera();
aCamera.Dolly(1.5);

// Set a background color for the renderer and set the size of the
// render window (expressed in pixels).
renWin.GetRenderer().SetBackground(1, 1, 1);
VtkPanelUtil.setSize(renWin, 640, 480);

// Note that when camera movement occurs (as it does in the Dolly()
// method), the clipping planes often need adjusting. Clipping planes
// consist of two planes: near and far along the view direction. The
// near plane clips out objects in front of the plane the far plane
// clips out objects behind the plane. This way only what is drawn
// between the planes is actually rendered.
renWin.GetRenderer().ResetCameraClippingRange();

// Setup panel
setLayout(new BorderLayout());
add(renWin, BorderLayout.CENTER);
}

public vtkPanel getRenWin() {
    return renWin;
}

```

```

    }

    public static void main(String s[]) {
        if (s.length == 0){
            return; //need a filename here
        }
        File theFile = new File(s[0]);
        //File theFile = new
            File("/Users/mmroden/Documents/MVSDownloadDirectory/Documents/1.2.840.113704.1.111.3384.1271766367.5/");
        AWTMedical3 panel = new AWTMedical3(theFile);

        JFrame frame = new JFrame("AWTMedical3");
        frame.setDefaultCloseOperation(JFrame.EXIT_ON_CLOSE);
        frame.getContentPane().add("Center", panel);
        frame.pack();
        frame.setVisible(true);
    }
}

```

## 29.2 BasicAnonymizer.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/BasicAnonymizer.exe gdcmData/012345.002.050.dcm out.dcm
 */
using System;
using gdcm;

public class MyWatcher : SimpleSubjectWatcher
{
    public MyWatcher(Subject s):base(s,"Override String"){
        protected override void StartFilter() {
            System.Console.WriteLine( "This is my start" );
        }
        protected override void EndFilter(){
            System.Console.WriteLine( "This is my end" );
        }
        protected override void ShowProgress(Subject caller, Event evt){
            ProgressEvent pe = ProgressEvent.Cast(evt);
            System.Console.WriteLine( "This is my progress: " + pe.GetProgress() );
        }
        protected override void ShowIteration(){
            System.Console.WriteLine( "This is my iteration" );
        }
        protected override void ShowAnonymization(Subject caller, Event evt){
/*
 * A couple of explanation are necessary here to understand how SWIG work
 * http://www.swig.org/Doc1.3/Java.html#adding_downcasts
 *
 * System.Console.WriteLine( "This is my Anonymization. Type: " + evt.GetEventName() );
 * System.Type type = evt.GetType();
 * System.Console.WriteLine( "This is my Anonymization. System.Type: " + type.ToString() );
 * System.Console.WriteLine( "This is my Anonymization. CheckEvent: " + ae.CheckEvent( evt ) );
 * System.Console.WriteLine( "This is my Anonymization. Processing Tag #" + ae.GetTag().toString() );
 */
            AnonymizeEvent ae = AnonymizeEvent.Cast(evt);
            if( ae != null )
            {
                Tag t = ae.GetTag();
                System.Console.WriteLine( "This is my Anonymization. Processing Tag #" + t.toString() );
            }
        }
    }
}

```

```

    }
    else
    {
        System.Console.WriteLine( "This is my Anonymization. Unhandled Event type: " + evt.GetEventName() );
    }
}
protected override void ShowAbort(){
    System.Console.WriteLine( "This is my abort" );
}
}

public class BasicAnonymizer
{
    public static int Main(string[] args)
    {
        gdcm.Global global = gdcm.Global.GetInstance();
        if( !global.LoadResourcesFiles() )
        {
            System.Console.WriteLine( "Could not LoadResourcesFiles" );
            return 1;
        }

        string file1 = args[0];
        string file2 = args[1];
        Reader reader = new Reader();
        reader.SetFileName( file1 );
        bool ret = reader.Read();
        if( !ret )
        {
            return 1;
        }

        string certpath = gdcm.Filename.Join(gdcm.Testing.
            GetSourceDirectory(), "/Testing/Source/Data/certificate.pem" );
        gdcm.CryptoFactory fact = gdcm.CryptoFactory.
            GetFactoryInstance();
        gdcm.CryptographicMessageSyntax cms = fact.CreateCMSProvider();
        if( !cms.ParseCertificateFile( certpath ) )
        {
            return 1;
        }

        //Anonymizer ano = new Anonymizer();
        SmartPtrAno sano = Anonymizer.New();
        Anonymizer ano = sano.__ref__();

        //SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(ano, "Anonymizer");
        MyWatcher watcher = new MyWatcher(ano);

        ano.SetFile( reader.GetFile() );
        ano.SetCryptographicMessageSyntax( cms );
        if( !ano.BasicApplicationLevelConfidentialityProfile() )
        {
            return 1;
        }

        Writer writer = new Writer();
        writer.SetFileName( file2 );
        writer.SetFile( ano.GetFile() );
        ret = writer.Write();
        if( !ret )
        {
            return 1;
        }

        return 0;
    }
}

```

## 29.3 BasicImageAnonymizer.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.

```

See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the above copyright notice for more information.

```

=====*/

/*
*/
using System;
using gdcm;

public class BasicImageAnonymizer
{
    public static int Main(string[] args)
    {
        string filename = args[0];

        // instanciate the reader:
        gdcm.ImageReader reader = new gdcm.ImageReader();
        reader.SetFileName( filename );

        if (!reader.Read()) return 1;

        Image ir = reader.GetImage();

        uint[] dims = {0, 0, 0};
        dims[0] = ir.GetDimension(0);
        dims[1] = ir.GetDimension(1);
        dims[2] = ir.GetDimension(2);
        System.Console.WriteLine( "Dim:" + dims[0] );
        System.Console.WriteLine( "Dim:" + dims[1] );
        System.Console.WriteLine( "Dim:" + dims[2] );

        // buffer to get the pixels
        byte[] buffer = new byte[ ir.GetBufferLength()];
        System.Console.WriteLine( "Dim:" + ir.GetBufferLength() );
        ir.GetBuffer( buffer );

        for (uint z = 0; z < dims[2]; z++)
        {
            for (uint y = 0; y < dims[1] / 2; y++) // only half Y
            {
                for (uint x = 0; x < dims[0] / 2; x++) // only half X
                {
                    buffer[ (z * dims[1] + y) * dims[0] + x ] = 0; // works when pixel type == UINT8
                }
            }
        }

        DataElement pixeldata = new DataElement( new Tag(0x7fe0,0x0010) );
        pixeldata.SetByteValue( buffer, new VL( (uint)buffer.Length ) );
        ir.SetDataElement( pixeldata );
        ir.SetTransferSyntax( new TransferSyntax( TransferSyntax.TSType.ExplicitVRLittleEndian ) );

        ImageChangeTransferSyntax change = new ImageChangeTransferSyntax();
        change.SetTransferSyntax( new TransferSyntax( TransferSyntax.TSType.JPEGLSLossless ) );
        change.SetInput( ir );
        if ( !change.Change() )
        {
            System.Console.WriteLine( "Could not change: " + filename );
            return 1;
        }

        ImageWriter writer = new ImageWriter();
        writer.SetFileName( "out.dcm" );
        writer.SetFile( reader.GetFile() );
        writer.SetImage( change.GetOutput() );
        bool ret = writer.Write();
        if ( !ret )
        {
            return 1;
        }

        return 0;
    }
}

```

## 29.4 CastConvertPhilips.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18     python --public /path/to/directory/
19 or
20     python --private /path/to/directory/
21
22     python --public --extension bak /path/to/directory/
23
24     rename -f 's/\.bak$//' *.bak
25
26 TODO:
27 http://docs.python.org/library/optparse.html#module-optparse
28 """
29
30 import vtkgdcm
31 import vtk
32 import sys
33 import gdcm
34
35 def ProcessOneFilePublic(filename, outfilename, tmpfile):
36     gdcm.ImageHelper.SetForceRescaleInterceptSlope(True)
37     vtkreader = vtkgdcm.vtkGDCMImageReader()
38     vtkreader.SetFileName( filename )
39     vtkreader.Update()
40
41     cast = vtk.vtkImageCast()
42     cast.SetInput( vtkreader.GetOutput() )
43     cast.SetOutputScalarTypeToUnsignedShort()
44
45     # vtkGDCMImageWriter does not support Sequence, so let's write a tmp file first:
46     # Some operation will actually be discarded (we simply need a temp storage)
47     vtkwriter = vtkgdcm.vtkGDCMImageWriter()
48     vtkwriter.SetFileName( tmpfile )
49     vtkwriter.SetMedicalImageProperties( vtkreader.GetMedicalImageProperties() )
50     vtkwriter.SetDirectionCosines( vtkreader.GetDirectionCosines() )
51     print "Format:", vtkreader.GetImageFormat()
52     vtkwriter.SetImageFormat( vtkreader.GetImageFormat() )
53     vtkwriter.SetInput( cast.GetOutput() )
54     #vtkwriter.Update()
55     vtkwriter.Write()
56
57     # ok now rewrite the exact same file as the original (keep all info)
58     # but use the Pixel Data Element from the written file
59     tmpreader = gdcm.ImageReader()
60     tmpreader.SetFileName( tmpfile )
61     if not tmpreader.Read():
62         sys.exit(1)
63
64     reader = gdcm.Reader()
65     reader.SetFileName( filename )
66     if not reader.Read():
67         sys.exit(1)
68
69     # Make sure to remove Slope/Rescale to avoid re-execution
70     ds = reader.GetFile().GetDataSet()
71     tags = [
72         gdcm.Tag(0x0028,0x1052),
73         gdcm.Tag(0x0028,0x1053),
74         gdcm.Tag(0x0028,0x1053),
75     ]
76     for tag in tags:
77         ds.Remove( tag )
78

```

```

79 writer = gdcmm.ImageWriter()
80 writer.SetFileName( outfile )
81 # Pass image from vtk written file
82 writer.SetImage( tmpreader.GetImage() )
83 # pass dataset from initial 'reader'
84 writer.SetFile( reader.GetFile() )
85 if not writer.Write():
86     sys.exit(1)
87
88 def ProcessOneFilePrivate(filename, outfile, tmpfile):
89     vtkreader = vtkgdcmm.vtkGDCMImageReader()
90     vtkreader.SetFileName( filename )
91     vtkreader.Update()
92
93
94 # (2005,1409)      DS      4      0.0
95 # (2005,140a)      DS      16     1.52283272283272
96
97 # (2005,0014)      LO      26     Philips MR Imaging DD 005
98 tag1 = gdcmm.PrivateTag(0x2005,0x09,"Philips MR Imaging DD 005")
99 tag2 = gdcmm.PrivateTag(0x2005,0x0a,"Philips MR Imaging DD 005")
100
101
102
103 # Need to access some private tags, reread the file (for now):
104 reader = gdcmm.Reader()
105 reader.SetFileName( filename )
106 if not reader.Read():
107     sys.exit(1)
108
109 ds = reader.GetFile().GetDataSet()
110
111 el1 = ds.GetDataElement( tag1 )
112 el2 = ds.GetDataElement( tag2 )
113
114
115 #pf = gdcmm.PythonFilter()
116 #pf.SetFile( reader.GetFile() )
117 #print el1.GetTag()
118
119 print el1.GetByteValue()
120 v1 = eval(el1.GetByteValue().GetBuffer())
121 print el2.GetByteValue()
122 v2 = eval(el2.GetByteValue().GetBuffer())
123
124 print v1
125 shift = v1
126 print v2
127 scale = v2
128
129 ss = vtk.vtkImageShiftScale()
130 ss.SetInput( vtkreader.GetOutput() )
131 # because VTK image shift / scale convention is inverted from DICOM make sure shift is 0
132 assert shift == 0
133 ss.SetShift( shift )
134 ss.SetScale( scale )
135 ss.SetOutputScalarTypeToUnsignedShort()
136 ss.Update()
137
138 # vtkGDCMImageWriter does not support Sequence, so let's write a tmp file first:
139 # Some operation will actually be discarded (we simply need a temp storage)
140 vtkwriter = vtkgdcmm.vtkGDCMImageWriter()
141 vtkwriter.SetFileName( tmpfile )
142 vtkwriter.SetMedicalImageProperties( vtkreader.GetMedicalImageProperties() )
143 vtkwriter.SetDirectionCosines( vtkreader.GetDirectionCosines() )
144 vtkwriter.SetImageFormat( reader.GetImageFormat() )
145 # do not pass shift/scale again
146 vtkwriter.SetInput( ss.GetOutput() )
147 #vtkwriter.Update()
148 vtkwriter.Write()
149
150 # ok now rewrite the exact same file as the original (keep all info)
151 # but use the Pixel Data Element from the written file
152 tmpreader = gdcmm.ImageReader()
153 tmpreader.SetFileName( tmpfile )
154 if not tmpreader.Read():
155     sys.exit(1)
156
157 writer = gdcmm.ImageWriter()
158 writer.SetFileName( outfile )
159 # Pass image from vtk written file

```

```

160 writer.SetImage( tmpreader.GetImage() )
161 # pass dataset from initial 'reader'
162 writer.SetFile( reader.GetFile() )
163 if not writer.Write():
164     sys.exit(1)
165
166 if __name__ == "__main__":
167
168     gdcmm.Trace.DebugOff()
169     gdcmm.Trace.WarningOff()
170     #filename = sys.argv[1]
171     #outfilename = sys.argv[2]
172     tmpfile = "/tmp/philips_rescaled.dcm"
173     #ProcessOneFile( filename, outfilename, tmpfile )
174     rescaletype = sys.argv[1]
175     assert rescaletype == "--public" or rescaletype == "--private"
176     dirname = sys.argv[2]
177     d = gdcmm.Directory()
178     d.Load( dirname )
179
180     for f in d.GetFilenames():
181         #print f
182         ProcessOneFilePublic( f, f + ".bak", tmpfile )
183
184
185 print "success"

```

## 29.5 ChangePrivateTags.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmmReader.h"
#include "gdcmmWriter.h"
#include "gdcmmPrivateTag.h"

int main(int argc, char* argv[] )
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " path/to/05148044-mr-siemens-avanto-syngo.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcmm::Reader reader;
    reader.SetFileName( filename );
    if (! reader.Read() )
    {
        return 1;
    }

    // (0029,0010) LO [SIEMENS CSA HEADER] # 18,1 Private Creator
    // (0029,0011) LO [SIEMENS MEDCOM HEADER ] # 22,1 Private Creator
    // (0029,0012) LO [SIEMENS MEDCOM HEADER2] # 22,1 Private Creator
    // [...]
    // (0029,1018) CS [MR] # 2,1 CSA Series Header Type
    // (0029,1134) CS [DB TO DICOM ] # 12,1 PMTF Information 4
    // (0029,1260) LO [com ] # 4,1 Series Workflow Status

    gdcmm::File &file = reader.GetFile();
    gdcmm::DataSet &ds = file.GetDataSet();

    // Declare private tag we need to find:
    gdcmm::PrivateTag pt1( 0x29,0x18, "SIEMENS CSA HEADER" );
    gdcmm::PrivateTag pt2( 0x29,0x34, "SIEMENS MEDCOM HEADER" );

```



```

gdcmm::PrivateTag pt3( 0x29,0x60, "SIEMENS MEDCOM HEADER2" );

const char str1[] = "GDCM was here 3!";
if( !ds.FindDataElement( pt1 ) ) return 1;
gdcmm::DataElement de1 = ds.GetDataElement( pt1 ); // Convert Private tag,
    into actual DataElement
std::cout << de1 << std::endl;
de1.SetByteValue( str1, (uint32_t)strlen(str1) );
ds.Replace( de1 );

const char str2[] = "GDCM was here 2!";
if( !ds.FindDataElement( pt2 ) ) return 1;
gdcmm::DataElement de2 = ds.GetDataElement( pt2 );
std::cout << de2 << std::endl;
de2.SetByteValue( str2, (uint32_t)strlen(str2) );
ds.Replace( de2 );

const char str3[] = "GDCM was here 3!";
if( !ds.FindDataElement( pt3 ) ) return 1;
gdcmm::DataElement de3 = ds.GetDataElement( pt3 );
std::cout << de3 << std::endl;
de3.SetByteValue( str3, (uint32_t)strlen(str3) );
ds.Replace( de3 );

gdcmm::Writer writer;
writer.SetFile( file );
writer.SetFileName( outfilename );
if ( !writer.Write() )
{
    return 1;
}

return 0;
}

```

## 29.6 ChangeSequenceUltrasound.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
#include "gdcmmReader.h"
#include "gdcmmWriter.h"
#include "gdcmmSmartPointer.h"
#include "gdcmmDataSetHelper.h"

/*
./ChangeSequenceUltrasound gdcmmData/D_CLUNIE_CT1_J2KI.dcm myoutput.dcm

This is the exact C++ translation of the original python example: ManipulateSequence.py
*/

int main(int argc, char* argv[] )
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcmm::Reader reader;
    reader.SetFileName( filename );
    if ( ! reader.Read() )
    {
        return 1;
    }
}

```

```

gdcM::File &file = reader.GetFile();
gdcM::DataSet &ds = file.GetDataSet();
gdcM::Tag tsis(0x0008,0x2112); // SourceImageSequence
if ( ds.FindDataElement( tsis ) )
{
    const gdcM::DataElement &sis = ds.GetDataElement( tsis );
    gdcM::SmartPointer<gdcM::SequenceOfItems> sqsis = sis.
        GetValueAsSQ();
    if ( sqsis && sqsis->GetNumberOfItems() )
    {
        gdcM::Item &item1 = sqsis->GetItem(1);
        gdcM::DataSet &nestedds = item1.GetNestedDataSet();
        gdcM::Tag tprcs(0x0040,0xa170); // PurposeOfReferenceCodeSequence
        if( nestedds.FindDataElement( tprcs ) )
        {
            const gdcM::DataElement &prcs = nestedds.GetDataElement( tprcs );
            gdcM::SmartPointer<gdcM::SequenceOfItems> sqprcs = prcs.
                GetValueAsSQ();
            if ( sqprcs && sqprcs->GetNumberOfItems() )
            {
                gdcM::Item &item2 = sqprcs->GetItem(1);
                gdcM::DataSet &nestedds2 = item2.GetNestedDataSet();
                // (0008,0104) LO [Uncompressed predecessor] # 24, 1 CodeMeaning
                gdcM::Tag tcm(0x0008,0x0104);
                if( nestedds2.FindDataElement( tcm ) )
                {
                    gdcM::DataElement cm = nestedds2.GetDataElement( tcm );
                    std::string mystr = "GDCM was here";
                    cm.SetByteValue( mystr.c_str(), (uint32_t)mystr.size() );
                    nestedds2.Replace( cm );
                }
            }
        }
    }
}

gdcM::Writer writer;
writer.SetFile( file );
writer.SetFileName( outfilename );
if ( !writer.Write() )
{
    return 1;
}

return 0;
}

```

## 29.7 CheckBigEndianBug.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * WARNING: This is a dev tool, do not use !
 *
 * Usage: after a gdcMconv, you would like to know if the conversion process is acceptable
 * sometime a vbindiff is acceptable, sometime it is not. In the case of the famous Philips
 * Little/Big Endian Explicit Transfer Syntax it is not easy to compare two files. However
 * this only impact byte ordering, thus we can compute byte-independant information to still
 * compare the files.
 */

#include "gdcMImageReader.h"
#include "gdcMImage.h"
#include "gdcMWriter.h"
#include "gdcMAttribute.h"

```

```

#include "gdcmSystem.h"

#include <iostream>
#include <fstream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input1.dcm input2.dcm" << std::endl;
        return 1;
    }
    const char *filename1 = argv[1];
    const char *filename2 = argv[2];

    gdcm::ImageReader reader1;
    reader1.SetFileName( filename1 );
    if( !reader1.Read() )
    {
        std::cerr << "Could not read: " << filename1 << std::endl;
        return 1;
    }

    gdcm::ImageReader reader2;
    reader2.SetFileName( filename2 );
    if( !reader2.Read() )
    {
        std::cerr << "Could not read: " << filename2 << std::endl;
        return 1;
    }

    // TODO: need a DataSet== operator implementation

    std::cout << "Both files can be read and looks like DICOM" << std::endl;

    size_t s1 = gdcm::System::FileSize(filename1);
    size_t s2 = gdcm::System::FileSize(filename2);

    if( s1 != s2 )
    {
        std::cout << "Size mismatch: " << s1 << " != " << s2 << std::endl;
        return 1;
    }
    else
    {
        std::cout << "Size match: " << s1 << " = " << s2 << std::endl;
    }

    std::ifstream is1( filename1, std::ios::binary );
    char *buffer1 = new char[s1];
    is1.read(buffer1, s1);

    std::ifstream is2( filename2, std::ios::binary );
    char *buffer2 = new char[s2];
    is2.read(buffer2, s2);

    assert( s1 == s2 );
    if( memcmp(buffer1, buffer2, s1 ) == 0 )
    {
        std::cout << "memcmp succeed ! File are bit identical" << std::endl;
    }
    else
    {
        std::cout << "memcmp failed!" << std::endl;
    }

    // Hum...memcmp failed, for big endian/ little endian inversion the histogram of bytes
    // should still be the same. So let's compute it
    // buffer2[0] = 1; // let's make the test fail
    std::multiset<char> set1( buffer1, buffer1 + s1 );
    std::multiset<char> set2( buffer2, buffer2 + s2 );

    if( set1 == set2 )
    {
        std::cout << "set1 == set2. Byte histogram seems valid" << std::endl;
    }
    else
    {
        std::cout << "set1 != set2" << std::endl;
    }
}

```

```

delete[] buffer1;
delete[] buffer2;

return 0;
}

```

## 29.8 ClinicalTrialAnnotate.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

/*
 * Dummy implementation of C.7.1.3 Clinical Trial Subject Module
 *
 * Usage:
 * ClinicalTrialAnnotate gdcmlData/012345.002.050.dcm out.dcm
 */

#include "gdcmlReader.h"
#include "gdcmlWriter.h"
#include "gdcmlAnonymizer.h"

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcml::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }

    // The output of gdcml::Reader is a gdcml::File
    //gdcml::File &file = reader.GetFile();

    // the dataset is the the set of element we are interested in:
    //gdcml::DataSet &ds = file.GetDataSet();

    gdcml::Anonymizer ano;
    ano.SetFile( reader.GetFile() );
    ano.RemoveGroupLength();
    ano.RemovePrivateTags();

    // PS 3.3 - 2008
    // C.7.1.3 Clinical Trial Subject Module
    // <entry group="0012" element="0010" vr="LO" vm="1" name="Clinical Trial Sponsor Name"/>
    ano.Replace( gdcml::Tag(0x12,0x10), "BigCompany name" );
    // <entry group="0012" element="0020" vr="LO" vm="1" name="Clinical Trial Protocol ID"/>
    ano.Replace( gdcml::Tag(0x12,0x20), "My Clinical Trial Protocol ID" );
    // <entry group="0012" element="0021" vr="LO" vm="1" name="Clinical Trial Protocol Name"/>
    ano.Replace( gdcml::Tag(0x12,0x21), "My Clinical Trial Protocol Name" );
    // <entry group="0012" element="0030" vr="LO" vm="1" name="Clinical Trial Site ID"/>
    ano.Replace( gdcml::Tag(0x12,0x30), "My Clinical Trial Site ID" );
    // <entry group="0012" element="0031" vr="LO" vm="1" name="Clinical Trial Site Name"/>
    ano.Replace( gdcml::Tag(0x12,0x31), "My Clinical Trial Site Name" );
    // <entry group="0012" element="0040" vr="LO" vm="1" name="Clinical Trial Subject ID"/>
    ano.Replace( gdcml::Tag(0x12,0x40), "My Clinical Trial Subject ID" );
    // <entry group="0012" element="0042" vr="LO" vm="1" name="Clinical Trial Subject Reading ID"/>

```

```

ano.Replace( gdcM:Tag(0x12,0x42), "My Clinical Trial Subject Reading ID" );

gdcM:Writer writer;
writer.SetFile( reader.GetFile() );
writer.SetFileName( outfilename );
if( !writer.Write() )
{
    return 1;
}

return 0;
}

```

## 29.9 ClinicalTrialIdentificationWorkflow.cs

This is a C# example on how to use Anonymizer

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * Typical usage on UNIX:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcM/debug-gcc/bin
 * $ mono bin/ClinicalTrialIdentificationWorkflow.exe input_dir output_dir
 */
using System;
using gdcM;

public class MyWatcher : SimpleSubjectWatcher
{
    public MyWatcher(Subject s):base(s,"Override String"){
        protected override void StartFilter() {
            System.Console.WriteLine( "This is my start" );
        }
        protected override void EndFilter(){
            System.Console.WriteLine( "This is my end" );
        }
        protected override void ShowProgress(Subject caller, Event evt){
            ProgressEvent pe = ProgressEvent.Cast(evt);
            System.Console.WriteLine( "This is my progress: " + pe.GetProgress() );
        }
        protected override void ShowIteration(){
            System.Console.WriteLine( "This is my iteration" );
        }
        protected override void ShowAnonymization(Subject caller, Event evt){
/*
 * A couple of explanation are necessary here to understand how SWIG work
 * http://www.swig.org/Doc1.3/Java.html#adding_downcasts
 *
 * System.Console.WriteLine( "This is my Anonymization. Type: " + evt.GetEventName() );
 * System.Type type = evt.GetType();
 * System.Console.WriteLine( "This is my Anonymization. System.Type: " + type.ToString() );
 * System.Console.WriteLine( "This is my Anonymization. CheckEvent: " + ae.CheckEvent( evt ) );
 * System.Console.WriteLine( "This is my Anonymization. Processing Tag #" + ae.GetTag().ToString() );
 */
            AnonymizeEvent ae = AnonymizeEvent.Cast(evt);
            if( ae != null )
            {
                Tag t = ae.GetTag();
                System.Console.WriteLine( "This is my Anonymization. Processing Tag #" + t.ToString() );
            }
            else
            {
                System.Console.WriteLine( "This is my Anonymization. Unhandled Event type: " + evt.GetEventName() );
            }
        }
    }
}

```

```

    }
    protected override void ShowAbort(){
        System.Console.WriteLine( "This is my abort" );
    }
}

public class ClinicalTrialIdentificationWorkflow
{
    public static bool ProcessOneFile( gdcm.Anonymizer ano , string filename, string
        outfilename )
    {
        Reader reader = new Reader();
        reader.SetFileName( filename );
        bool ret = reader.Read();
        if( !ret )
        {
            return false;
        }
        // Pass in the file:
        ano.SetFile( reader.GetFile() );

        // First step, let's protect all Patient information as per
        // PS 3.15 / E.1 / Basic Application Level Confidentiality Profile
        if( !ano.BasicApplicationLevelConfidentialityProfile() )
        {
            return false;
        }

        // Now let's pass in all Clinical Trial fields
        // PS 3.3 - 2008 / C.7.1.3 Clinical Trial Subject Module
        /*
        Clinical Trial Sponsor Name (0012,0010) 1 The name of the clinical trial sponsor. See C.7.1.3.1.1.
        Clinical Trial Protocol ID (0012,0020) 1 Identifier for the noted protocol. See C.7.1.3.1.2.
        Clinical Trial Protocol Name (0012,0021) 2 The name of the clinical trial protocol. See C.7.1.3.1.3.
        Clinical Trial Site ID (0012,0030) 2 The identifier of the site responsible for submitting clinical
            trial data. See C.7.1.3.1.4.
        Clinical Trial Site Name (0012,0031) 2 Name of the site responsible for submitting clinical trial data.
            See C.7.1.3.1.5
        Clinical Trial Subject ID (0012,0040) 1C The assigned identifier for the clinical trial subject. See
            C.7.1.3.1.6. Shall be present if Clinical Trial Subject Reading ID (0012,0042) is absent. May be present
            otherwise.
        Clinical Trial Subject Reading ID (0012,0042) 1C Identifies the subject for blinded evaluations. Shall
            be present if Clinical Trial Subject ID (0012,0040) is absent. May be present otherwise. See C.7.1.3.1.7.
        */
        ano.Replace( new gdcm.Tag(0x0012,0x0010), "MySponsorName");
        ano.Replace( new gdcm.Tag(0x0012,0x0020), "MyProtocolID");
        ano.Replace( new gdcm.Tag(0x0012,0x0021), "MyProtocolName");
        ano.Replace( new gdcm.Tag(0x0012,0x0030), "MySiteId");
        ano.Replace( new gdcm.Tag(0x0012,0x0031), "MySiteName");
        ano.Replace( new gdcm.Tag(0x0012,0x0040), "MySponsorId");
        ano.Replace( new gdcm.Tag(0x0012,0x0050), "MyTPId");
        ano.Replace( new gdcm.Tag(0x0012,0x0051), "MyTPDescription");

        // The following two are not required as they are guaranteed to be filled in by the
        // Basic Application Level Confidentiality Profile. Only override if you understand what
        // you are doing
        //ano.Replace( new gdcm.Tag(0x0012,0x0062), "YES");
        //ano.Replace( new gdcm.Tag(0x0012,0x0063), "My Super Duper Anonymization Overload");

        // We might be generating a subdirectory. Let's make sure the subdir exist:
        gdcm.FileName fn = new gdcm.FileName( outfilename );
        string subdir = fn.GetPath();
        if( !gdcm.PosixEmulation.MakeDirectory( subdir ) )
        {
            return false;
        }

        gdcm.FileMetaInformation fmi = ano.GetFile().GetHeader();
        // The following three lines make sure to regenerate any value:
        fmi.Remove( new gdcm.Tag(0x0002,0x0012) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0013) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0016) );

        Writer writer = new Writer();
        writer.SetFileName( outfilename );
        writer.SetFile( ano.GetFile() );
        ret = writer.Write();
        if( !ret )
        {
            return false;
        }
    }
}

```

```

    return true;
}

public static int Main(string[] args)
{
    gdcmm.FileMetaInformation.
        SetSourceApplicationEntityTitle( "My ClinicalTrial App" );

    // http://www.oid-info.com/get/1.3.6.1.4.17434
    string THERALYS_ORG_ROOT = "1.3.6.1.4.17434";
    gdcmm.UIDGenerator.SetRoot( THERALYS_ORG_ROOT );
    System.Console.WriteLine( "Root dir is now: " + gdcmm.UIDGenerator.
        GetRoot() );

    gdcmm.Global global = gdcmm.Global.GetInstance();
    if( !global.LoadResourcesFiles() )
    {
        System.Console.WriteLine( "Could not LoadResourcesFiles" );
        return 1;
    }

    if( args.Length != 2 )
    {
        System.Console.WriteLine( "Usage:" );
        System.Console.WriteLine( "ClinicalTrialIdentificationWorkflow input_dir output_dir" );
        return 1;
    }

    string dir1 = args[0];
    string dir2 = args[1];

    // Check input is valid:
    if( !gdcmm.PosixEmulation.FileIsDirectory(dir1) )
    {
        System.Console.WriteLine( "Input directory: " + dir1 + " does not exist. Sorry" );
        return 1;
    }
    if( !gdcmm.PosixEmulation.FileIsDirectory(dir2) )
    {
        System.Console.WriteLine( "Output directory: " + dir2 + " does not exist. Sorry" );
        return 1;
    }

    // Recursively search all file within this toplevel directory:
    Directory d = new Directory();
    uint nfiles = d.Load( dir1, true );
    if(nfiles == 0) return 1;

    // Let's use the pre-shipped certificate of GDCM.
    string certpath = gdcmm.FileName.Join(gdcmm.Testing.
        GetSourceDirectory(), "/Testing/Source/Data/certificate.pem" );
    gdcmm.CryptoFactory fact = gdcmm.CryptoFactory.
        GetFactoryInstance();
    gdcmm.CryptographicMessageSyntax cms = fact.CreateCMSProvider();
    if( !cms.ParseCertificateFile( certpath ) )
    {
        System.Console.WriteLine( "PEM Certificate : " + certpath + " could not be read. Sorry" );
        return 1;
    }

    //Anonymizer ano = new Anonymizer();
    // A reference to an actual C++ instance is required here:
    SmartPtrAno sano = Anonymizer.New();
    Anonymizer ano = sano.__ref__();

    //SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(ano, "Anonymizer");
    MyWatcher watcher = new MyWatcher(ano);

    // Explicitely specify the Cryptographic Message Syntax to use:
    ano.SetCryptographicMessageSyntax( cms );

    // Process all filenames:
    FilenamesType filenames = d.GetFilesNames();
    for( uint i = 0; i < nfiles; ++i )
    {
        string filename = filenames[ (int)i ];
        string outfilename = filename.Replace( dir1, dir2 );
        System.Console.WriteLine( "Filename: " + filename );
        System.Console.WriteLine( "Out Filename: " + outfilename );
        if( !ProcessOneFile( ano , filename, outfilename ) )
        {

```

```

        System.Console.WriteLine( "Could not process filename: " + filename );
        return 1;
    }
}

return 0;
}
}

```

## 29.10 CompressImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 *
 */

#include "gdcmImageReader.h"
#include "gdcmImage.h"
#include "gdcmWriter.h"
#include "gdcmAttribute.h"
#include "gdcmImageWriter.h"
#include "gdcmImageChangeTransferSyntax.h"

#include <iostream>
#include <fstream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcm::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }

    // The output of gdcm::Reader is a gdcm::File
    //gdcm::File &file = reader.GetFile();

    // the dataset is the the set of element we are interested in:
    //gdcm::DataSet &ds = file.GetDataSet();

    const gdcm::Image &image = reader.GetImage();
    image.Print( std::cout );

    gdcm::ImageChangeTransferSyntax change;
    change.SetTransferSyntax(
        gdcm::TransferSyntax::JPEG2000Lossless );
    change.SetTransferSyntax(
        gdcm::TransferSyntax::JPEGLosslessProcess14_1 );
    //change.SetTransferSyntax( gdcm::TransferSyntax::JPEGBaselineProcess1 );
    //change.SetTransferSyntax( image.GetTransferSyntax() );
    change.SetInput( image );
    bool b = change.Change();
    if( !b )
    {
        std::cerr << "Could not change the Transfer Syntax" << std::endl;
    }
}

```



```

        return 1;
    }

    //std::ofstream out( outfilename, std::ios::binary );
    //image.GetBuffer2(out);
    //out.close();
    gdcm::ImageWriter writer;
    writer.SetImage( change.GetOutput() );
    writer.SetFile( reader.GetFile() );
    writer.SetFileName( outfilename );
    if( !writer.Write() )
    {
        return 1;
    }

    return 0;
}

```

## 29.11 CompressLossyJPEG.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Perso/gdcm/debug-gcc/bin
 * $ mono bin/CompressLossyJPEG.exe input.dcm output.dcm
 */

using System;
using gdcm;

public class CompressLossyJPEG
{
    public static int Main(string[] args)
    {
        if( args.Length < 2 )
        {
            System.Console.WriteLine( " input.dcm output.dcm" );
            return 1;
        }
        string filename = args[0];
        string outfilename = args[1];

        ImageReader reader = new ImageReader();
        reader.SetFileName( filename );
        if( !reader.Read() )
        {
            System.Console.WriteLine( "Could not read: " + filename );
            return 1;
        }

        // The output of gdcm::Reader is a gdcm::File
        File file = reader.GetFile();

        // the dataset is the the set of element we are interested in:
        DataSet ds = file.GetDataSet();

        Image image = reader.GetImage();
        //image.Print( cout );

        ImageChangeTransferSyntax change = new ImageChangeTransferSyntax();
        TransferSyntax targetts = new TransferSyntax( TransferSyntax.TSType.JPEGBaselineProcess1 );
        change.SetTransferSyntax( targetts );

        // Setup our JPEGCodec, warning it should be compatible with JPEGBaselineProcess1
        JPEGCodec jpegcodec = new JPEGCodec();

```

```

    if( !jpegcodec.CanCode( targetts ) )
    {
        System.Console.WriteLine( "Something went really wrong, JPEGCodec cannot handle JPEGBaselineProcess1"
        );
        return 1;
    }
    jpegcodec.SetLossless( false );
    jpegcodec.SetQuality( 50 ); // poor quality !
    change.SetUserCodec( jpegcodec ); // specify the codec to use to the ImageChangeTransferSyntax

    change.SetInput( image );
    bool b = change.Change();
    if( !b )
    {
        System.Console.WriteLine( "Could not change the Transfer Syntax" );
        return 1;
    }

    ImageWriter writer = new ImageWriter();
    writer.SetImage( (gdcm.Image)change.GetOutput() );
    writer.SetFile( reader.GetFile() );
    writer.SetFileName( outfilename );
    if( !writer.Write() )
    {
        System.Console.WriteLine( "Could not write: " + outfilename );
        return 1;
    }

    return 0;
}
}

```

## 29.12 Compute3DSpacing.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader2.h"
#include "vtkImageChangeInformation.h"
#include "vtkStringArray.h"
#include "gdcmIPPSorter.h"

/*
 * Simple example to check computation of spacing within vtkGDCMImageReader2
 * This is a direct implementation of:
 *
 * http://gdcm.sourceforge.net/wiki/index.php/
 *   Using_GDCM_API#Automatic_ordering_of_slices_for_vtkGDCMImageReader.SetFileNames
 *
 * For more advanced information on how 3D spacing is being computed see:
 *
 * - http://gdcm.sourceforge.net/html/classgdcm_1_1IPPSorter.html
 *
 * Usage:
 *
 * $ Compute3DSpacing SIEMENS_MAGNETOM-12-MONO2-FileSeq0.dcm \
 *   SIEMENS_MAGNETOM-12-MONO2-FileSeq1.dcm \
 *   SIEMENS_MAGNETOM-12-MONO2-FileSeq2.dcm \
 *   SIEMENS_MAGNETOM-12-MONO2-FileSeq3.dcm
 */

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;

    std::vector<std::string> filenames;

```

```

for( int i = 1; i < argc; ++i )
{
    filenames.push_back( argv[i] );
}

gdcmm::IPPSorter s;
s.SetComputeZSpacing( true );
s.SetZSpacingTolerance( 1e-3 );
bool b = s.Sort( filenames );
if( !b )
{
    std::cerr << "Failed to sort files" << std::endl;
    return 1;
}
std::cout << "Sorting succeeded:" << std::endl;
//s.Print( std::cout );

std::cout << "Found z-spacing:" << std::endl;
std::cout << s.GetZSpacing() << std::endl;
const double ippszspacing = s.GetZSpacing();

const std::vector<std::string> & sorted = s.GetFilenames();
vtkGDCMImageReader2 * reader = vtkGDCMImageReader2::New();
vtkStringArray *files = vtkStringArray::New();
std::vector< std::string >::const_iterator it = sorted.begin();
for( ; it != sorted.end(); ++it )
{
    const std::string &f = *it;
    files->InsertNextValue( f.c_str() );
}
reader->SetFileNames( files );
reader->Update();

const vtkFloatingPointType *spacing = reader->GetOutput()->GetSpacing();
vtkImageChangeInformation *v16 = vtkImageChangeInformation::New();
#if (VTK_MAJOR_VERSION >= 6)
v16->SetInputConnection( reader->GetOutputPort() );
#else
v16->SetInput( reader->GetOutput() );
#endif
v16->SetOutputSpacing( spacing[0], spacing[1], ippszspacing );
v16->Update();

v16->GetOutput()->Print( std::cout );

return 0;
}

```

## 29.13 Convert16BitsTo8Bits.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkImageData.h"
#include "vtkImageCast.h"

#include "gdcmmTesting.h"
// The following file is 16/16/15 but the scalar range of the image is [0,192]
// it could be safely stored as 8bits instead:
// gdcmmData/012345.002.050.dcm

int main(int, char *[])
{
    const char *directory = gdcmm::Testing::GetDataRoot();
    if(!directory) return 1;
}

```

```

std::string file = std::string(directory) + "/012345.002.050.dcm";
std::cout << file << std::endl;

vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
reader->SetFileName( file.c_str() );
reader->Update();
//reader->GetOutput()->Print( std::cout );

vtkImageCast *cast = vtkImageCast::New();
#if (VTK_MAJOR_VERSION >= 6)
    cast->SetInputConnection( reader->GetOutputPort() );
#else
    cast->SetInput( reader->GetOutput() );
#endif
cast->SetOutputScalarTypeToUnsignedChar();

vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
writer->SetFileName( "/tmp/cast.dcm" );
#if (VTK_MAJOR_VERSION >= 6)
    writer->SetInputConnection( cast->GetOutputPort() );
#else
    writer->SetInput( cast->GetOutput() );
#endif
writer->SetImageFormat( reader->GetImageFormat() );
writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
writer->SetDirectionCosines( reader->GetDirectionCosines() );
writer->SetShift( reader->GetShift() );
writer->SetScale( reader->GetScale() );
writer->Write();

reader->Delete();
cast->Delete();
writer->Delete();

return 0;
}

```

## 29.14 ConvertMPL.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 display a DICOM image with matplotlib via numpy
17
18 Caveats:
19 - Does not support UINT12/INT12
20
21 Usage:
22
23   python ConvertNumpy.py "IM000000"
24
25 Thanks:
26   plotting example - Ray Schumacher 2009
27 """
28
29 import gdcm
30 import numpy
31 from pylab import *
32
33
34 def get_gdcm_to_numpy_typemap():
35     """Returns the GDCM Pixel Format to numpy array type mapping."""
36     _gdcm_np = {gdcm.PixelFormat.UINT8 :numpy.int8,
37                 gdcm.PixelFormat.INT8  :numpy.uint8,

```

```

38         gdcM.PixelFormat.UINT16 :numpy.uint16,
39         gdcM.PixelFormat.INT16  :numpy.int16,
40         gdcM.PixelFormat.UINT32 :numpy.uint32,
41         gdcM.PixelFormat.INT32  :numpy.int32,
42         gdcM.PixelFormat.FLOAT32:numpy.float32,
43         gdcM.PixelFormat.FLOAT64:numpy.float64 }
44     return _gdcM_np
45
46 def get_numpy_array_type(gdcM_pixel_format):
47     """Returns a numpy array typecode given a GDCM Pixel Format."""
48     return get_gdcM_to_numpy_typemap()[gdcM_pixel_format]
49
50 def gdcM_to_numpy(image):
51     """Converts a GDCM image to a numpy array.
52     """
53     pf = image.GetPixelFormat().GetScalarType()
54     print 'pf', pf
55     print image.GetPixelFormat().GetScalarTypeAsString()
56     assert pf in get_gdcM_to_numpy_typemap().keys(), \
57         "Unsupported array type %s"%pf
58     d = image.GetDimension(0), image.GetDimension(1)
59     print 'Image Size: %d x %d' % (d[0], d[1])
60     dtype = get_numpy_array_type(pf)
61     gdcM_array = image.GetBuffer()
62     ## use float for accurate scaling
63     result = numpy.frombuffer(gdcM_array, dtype=dtype).astype(float)
64     ## optional gamma scaling
65     #maxV = float(result[result.argmax()])
66     #result = result + .5*(maxV-result)
67     #result = numpy.log(result+50) ## apprx background level
68     result.shape = d
69     return result
70
71 if __name__ == "__main__":
72     import sys
73     r = gdcM.ImageReader()
74     filename = sys.argv[1]
75     r.SetFileName( filename )
76     if not r.Read(): sys.exit(1)
77     numpy_array = gdcM_to_numpy( r.GetImage() )
78
79     subplot(111)# one plot, on left
80     title(filename)
81     ## many colormaps are available
82     imshow(numpy_array, interpolation='bilinear', cmap=cm.jet)
83     ## set the plot sizes and placement
84     subplots_adjust(bottom=0.1, right=0.8, top=0.9)
85     cax = axes([0.85, 0.1, 0.075, 0.8])
86     colorbar(cax=cax)
87     title('values')
88     get_current_fig_manager().window.title('plot')
89     show()

```

## 29.15 ConvertMultiFrameToSingleFrame.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkImageData.h"
#include "vtkStringArray.h"

#include "gdcMTesting.h"
#include "gdcMFilenameGenerator.h"

int main(int argc, char *argv[])

```

```

{
    std::string filename;
    if( argc <= 1 )
    {
        const char *directory = gdcm::Testing::GetDataRoot();
        if(!directory) return 1;
        std::string file = std::string(directory) + "/US-PAL-8-10x-echo.dcm";
        filename = file;
    }
    else
    {
        filename = argv[1];
    }
    std::cout << "file: " << filename << std::endl;

    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( filename.c_str() );
    reader->Update();
    //reader->GetOutput()->Print( std::cout );

    int dims[3];
    reader->GetOutput()->GetDimensions( dims );

    std::ostringstream os;
    os << "singleframe";
    os << "%04d.dcm";
    gdcm::FilenameGenerator fg;
    fg.SetPattern( os.str().c_str() );
    unsigned int nfiles = dims[2];
    fg.SetNumberOfFileNames( nfiles );
    bool b = fg.Generate();
    if( !b )
    {
        std::cerr << "FilenameGenerator::Generate() failed" << std::endl;
        return 1;
    }
    if( !fg.GetNumberOfFileNames() )
    {
        std::cerr << "FilenameGenerator::Generate() failed somehow..." << std::endl;
        return 1;
    }

    // By default write them as Secondary Capture (for portability)
    vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
    vtkStringArray *filenames = vtkStringArray::New();
    for(unsigned int i = 0; i < fg.GetNumberOfFileNames(); ++i)
    {
        filenames->InsertNextValue( fg.GetFilename(i) );
    }
    assert( filenames->GetNumberOfValues() == (int)fg.GetNumberOfFileNames() );
    writer->SetFileNames( filenames );
    filenames->Delete();
    writer->SetFileDimensionality( 2 );
    #if (VTK_MAJOR_VERSION >= 6)
    writer->SetInputConnection( reader->GetOutputPort() );
    #else
    writer->SetInput( reader->GetOutput() );
    #endif
    writer->SetImageFormat( reader->GetImageFormat() );
    writer->Write();

    reader->Delete();
    writer->Delete();

    return 0;
}

```

## 29.16 ConvertNumpy.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #

```

```

9 #      This software is distributed WITHOUT ANY WARRANTY; without even
10 #      the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #      PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 This module add support for converting a gdcm.Image to a numpy array.
17
18 Caveats:
19 - Does not support UINT12/INT12
20
21 Removed:
22 - float16 is defined in GDCM API but no implementation exist for it ...
23 """
24
25 import gdcm
26 import numpy
27
28 def get_gdcm_to_numpy_typemap():
29     """Returns the GDCM Pixel Format to numpy array type mapping."""
30     _gdcm_np = {gdcm.PixelFormat.UINT8 :numpy.int8,
31                 gdcm.PixelFormat.INT8 :numpy.uint8,
32                 #gdcm.PixelFormat.UINT12 :numpy.uint12,
33                 #gdcm.PixelFormat.INT12 :numpy.int12,
34                 gdcm.PixelFormat.UINT16 :numpy.uint16,
35                 gdcm.PixelFormat.INT16 :numpy.int16,
36                 gdcm.PixelFormat.UINT32 :numpy.uint32,
37                 gdcm.PixelFormat.INT32 :numpy.int32,
38                 #gdcm.PixelFormat.FLOAT16:numpy.float16,
39                 gdcm.PixelFormat.FLOAT32:numpy.float32,
40                 gdcm.PixelFormat.FLOAT64:numpy.float64 }
41     return _gdcm_np
42
43 def get_numpy_array_type(gdcm_pixel_format):
44     """Returns a numpy array typecode given a GDCM Pixel Format."""
45     return get_gdcm_to_numpy_typemap()[gdcm_pixel_format]
46
47 def gdcm_to_numpy(image):
48     """Converts a GDCM image to a numpy array.
49     """
50     pf = image.GetPixelFormat()
51
52     assert pf.GetScalarType() in get_gdcm_to_numpy_typemap().keys(), \
53         "Unsupported array type %s"%pf
54
55     shape = image.GetDimension(0) * image.GetDimension(1), pf.GetSamplesPerPixel()
56     if image.GetNumberOfDimensions() == 3:
57         shape = shape[0] * image.GetDimension(2), shape[1]
58
59     dtype = get_numpy_array_type(pf.GetScalarType())
60     gdcm_array = image.GetBuffer()
61     result = numpy.frombuffer(gdcm_array, dtype=dtype)
62     result.shape = shape
63     return result
64
65 if __name__ == "__main__":
66     import sys
67     r = gdcm.ImageReader()
68     filename = sys.argv[1]
69     r.SetFileName( filename )
70     if not r.Read():
71         sys.exit(1)
72
73     numpy_array = gdcm_to_numpy( r.GetImage() )
74     print numpy_array

```

## 29.17 ConvertPIL.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #

```

```

9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 save a DICOM image with PIL via numpy
17
18 Caveats:
19 - Does not support UINT12/INT12
20
21 Usage:
22
23 python ConvertNumpy.py "IM000000"
24
25 Thanks:
26 plotting example - Ray Schumacher 2009
27 """
28
29 import gdcm
30 import numpy
31 from PIL import Image, ImageOps
32
33
34 def get_gdcm_to_numpy_typemap():
35     """Returns the GDCM Pixel Format to numpy array type mapping."""
36     _gdcm_np = {gdcm.PixelFormat.UINT8 :numpy.int8,
37                 gdcm.PixelFormat.INT8 :numpy.uint8,
38                 gdcm.PixelFormat.UINT16 :numpy.uint16,
39                 gdcm.PixelFormat.INT16 :numpy.int16,
40                 gdcm.PixelFormat.UINT32 :numpy.uint32,
41                 gdcm.PixelFormat.INT32 :numpy.int32,
42                 gdcm.PixelFormat.FLOAT32 :numpy.float32,
43                 gdcm.PixelFormat.FLOAT64 :numpy.float64 }
44     return _gdcm_np
45
46 def get_numpy_array_type(gdcm_pixel_format):
47     """Returns a numpy array typecode given a GDCM Pixel Format."""
48     return get_gdcm_to_numpy_typemap()[gdcm_pixel_format]
49
50 def gdcm_to_numpy(image):
51     """Converts a GDCM image to a numpy array.
52     """
53     pf = image.GetPixelFormat().GetScalarType()
54     print 'pf', pf
55     print image.GetPixelFormat().GetScalarTypeAsString()
56     assert pf in get_gdcm_to_numpy_typemap().keys(), \
57         "Unsupported array type %s"%pf
58     d = image.GetDimension(0), image.GetDimension(1)
59     print 'Image Size: %d x %d' % (d[0], d[1])
60     dtype = get_numpy_array_type(pf)
61     gdcm_array = image.GetBuffer()
62     result = numpy.frombuffer(gdcm_array, dtype=dtype)
63     maxV = float(result[result.argmax()])
64     ## linear gamma adjust
65     #result = result + .5*(maxV-result)
66     ## log gamma
67     result = numpy.log(result+50) ## 50 is apprxx background level
68     maxV = float(result[result.argmax()])
69     result = result*(2.**8/maxV) ## histogram stretch
70     result.shape = d
71     return result
72
73 if __name__ == "__main__":
74     import sys
75     r = gdcm.ImageReader()
76     filename = sys.argv[1]
77     r.SetFileName( filename )
78     if not r.Read(): sys.exit(1)
79     numpy_array = gdcm_to_numpy( r.GetImage() )
80     ## L is 8 bit grey
81     ## http://www.pythonware.com/library/pil/handbook/concepts.htm
82     pilImage = Image.frombuffer('L',
83                                numpy_array.shape,
84                                numpy_array.astype(numpy.uint8),
85                                'raw','L',0,1)
86     ## cutoff removes background noise and spikes
87     pilImage = ImageOps.autocontrast(pilImage, cutoff=.1)
88     pilImage.save(sys.argv[1]+' .jpg')

```



## 29.18 ConvertRGBToLuminance.cxx

```

/*=====

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkImageData.h"
#include "vtkImageLuminance.h"

#include "gdcmTesting.h"

// There is no such thing as MR Image Storage + Photometric Interpretation = RGB
// let's rewrite that into a proper single component image:
int main(int, char *[])
{
    const char *directory = gdcm::Testing::GetDataRoot();
    if(!directory) return 1;
    std::string file = std::string(directory) + "/SIEMENS-MR-RGB-16Bits.dcm";
    std::cout << file << std::endl;

    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( file.c_str() );
    reader->Update();
    //reader->GetOutput()->Print( std::cout );

    vtkImageLuminance *luminance = vtkImageLuminance::New();
    #if (VTK_MAJOR_VERSION >= 6)
        luminance->SetInputConnection( reader->GetOutputPort() );
    #else
        luminance->SetInput( reader->GetOutput() );
    #endif

    vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
    writer->SetFileName( "/tmp/bla.dcm" );
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputConnection( luminance->GetOutputPort() );
    #else
        writer->SetInput( luminance->GetOutput() );
    #endif
    //writer->SetImageFormat( reader->GetImageFormat() ); // Do NOT pass image format
    writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
    writer->SetDirectionCosines( reader->GetDirectionCosines() );
    writer->SetShift( reader->GetShift() );
    writer->SetScale( reader->GetScale() );
    writer->Write();

    // TODO:
    //vtkImageAppendComponents.h

    reader->Delete();
    luminance->Delete();
    writer->Delete();

    return 0;
}

```

## 29.19 ConvertSingleBitTo8Bits.cxx

```

/*=====

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.

```

See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the above copyright notice for more information.

```

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkImageData.h"
#include "vtkImageCast.h"
#include "vtkPointData.h"
#include "vtkBitArray.h"
#include "vtkUnsignedCharArray.h"

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( filename );
    reader->Update();
    //reader->GetOutput()->Print( std::cout );

    vtkDataArray* array = reader->GetOutput()->GetPointData()->GetScalars();
    vtkBitArray *barray = vtkBitArray::SafeDownCast( array );
    if( !barray ) return false;
    vtkIdType nvalues = array->GetNumberOfTuples();
    vtkUnsignedCharArray *uarray = vtkUnsignedCharArray::New();
    uarray->SetNumberOfTuples( nvalues );
    for(vtkIdType i = 0; i < nvalues; ++i)
    {
        uarray->SetValue( i, (unsigned char)barray->GetValue(i) );
    }

    vtkImageData *copy = vtkImageData::New();
    //
    http://www.vtk.org/Wiki/VTK/VTK_6_Migration/Changes_to_Scalars_Manipulation_Functions#AllocateScalars.28.29
    copy->SetExtent( reader->GetOutput()->GetExtent() );
    #if (VTK_MAJOR_VERSION >= 6)
    copy->AllocateScalars(VTK_UNSIGNED_CHAR, 3);
    #else
    copy->SetScalarType( VTK_UNSIGNED_CHAR );
    copy->AllocateScalars();
    #endif

    //uarray->Print( std::cout );
    //copy->GetPointData()->GetScalars()->Print( std::cout );
    copy->GetPointData()->SetScalars( uarray );
    uarray->Delete();

    vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
    writer->SetFileName( outfile );
    //writer->SetInput( cast->GetOutput() );
    #if (VTK_MAJOR_VERSION >= 6)
    writer->SetInputData( copy );
    #else
    writer->SetInput( copy );
    #endif
    writer->SetImageFormat( reader->GetImageFormat() );
    writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
    writer->SetDirectionCosines( reader->GetDirectionCosines() );
    writer->SetShift( reader->GetShift() );
    writer->SetScale( reader->GetScale() );
    writer->SetFileDimensionality( reader->GetFileDimensionality() );
    writer->Write();

    reader->Delete();
    copy->Delete();
    writer->Delete();

    return 0;
}

```

## 29.20 ConvertToQImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example shows how to setup the pipeline from a gdcm::ImageReader into a
 * Qt QImage data structure.
 * It only handles 2D image.
 *
 * Ref:
 * http://doc.trolltech.com/4.5/qimage.html
 *
 * Usage:
 * ConvertToQImage gdcmData/012345.002.050.dcm output.png
 *
 * Thanks:
 * Sylvain ADAM (sylvain51 hotmail com) for contributing this example
 */

#include "gdcmImageReader.h"
#include <QImage>
#include <QImageWriter>

bool ConvertToFormat_RGB888(gdcm::Image const & gimage, char *buffer, QImage* &imageQt)
{
    const unsigned int* dimension = gimage.GetDimensions();

    unsigned int dimX = dimension[0];
    unsigned int dimY = dimension[1];

    gimage.GetBuffer(buffer);

    // Let's start with the easy case:
    if( gimage.GetPhotometricInterpretation() ==
        gdcm::PhotometricInterpretation::RGB )
    {
        if( gimage.GetPixelFormat() != gdcm::PixelFormat::UINT8 )
        {
            return false;
        }
        unsigned char *ubuffer = (unsigned char*)buffer;
        // QImage::Format_RGB888 13 The image is stored using a 24-bit RGB format (8-8-8).
        imageQt = new QImage((unsigned char *)ubuffer, dimX, dimY, 3*dimX, QImage::Format_RGB888);
    }
    else if( gimage.GetPhotometricInterpretation() ==
        gdcm::PhotometricInterpretation::MONOCHROME2 )
    {
        if( gimage.GetPixelFormat() == gdcm::PixelFormat::UINT8 )
        {
            // We need to copy each individual 8bits into R / G and B:
            unsigned char *ubuffer = new unsigned char[dimX*dimY*3];
            unsigned char *pubuffer = ubuffer;
            for(unsigned int i = 0; i < dimX*dimY; i++)
            {
                *pubuffer++ = *buffer;
                *pubuffer++ = *buffer;
                *pubuffer++ = *buffer;
            }

            imageQt = new QImage(ubuffer, dimX, dimY, QImage::Format_RGB888);
        }
        else if( gimage.GetPixelFormat() == gdcm::PixelFormat::INT16 )
        {
            // We need to copy each individual 16bits into R / G and B (truncate value)
            short *buffer16 = (short*)buffer;
            unsigned char *ubuffer = new unsigned char[dimX*dimY*3];
            unsigned char *pubuffer = ubuffer;
            for(unsigned int i = 0; i < dimX*dimY; i++)
            {

```

```

        // Scalar Range of gdcmlData/012345.002.050.dcm is [0,192], we could simply do:
        // *pubbuffer++ = *buffer16;
        // *pubbuffer++ = *buffer16;
        // *pubbuffer++ = *buffer16;
        // instead do it right:
        *pubbuffer++ = (unsigned char)std::min(255, (32768 + *buffer16) / 255);
        *pubbuffer++ = (unsigned char)std::min(255, (32768 + *buffer16) / 255);
        *pubbuffer++ = (unsigned char)std::min(255, (32768 + *buffer16) / 255);
        buffer16++;
    }

    QImage *imageQt = new QImage(ubuffer, dimX, dimY, QImage::Format_RGB888);
}
else
{
    std::cerr << "Pixel Format is: " << gimage.GetPixelFormat() << std::endl;
    return false;
}
}
else
{
    std::cerr << "Unhandled PhotometricInterpretation: " << gimage.
        GetPhotometricInterpretation() << std::endl;
    return false;
}

return true;
}

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcm::ImageReader ir;
    ir.SetFileName( filename );
    if(!ir.Read())
    {
        //Read failed
        return 1;
    }

    std::cout<<"Getting image from ImageReader..."<<std::endl;

    const gdcm::Image &gimage = ir.GetImage();
    std::vector<char> vbuffer;
    vbuffer.resize( gimage.GetBufferLength() );
    char *buffer = &vbuffer[0];

    QImage *imageQt = NULL;
    if( !ConvertToFormat_RGB888( gimage, buffer, imageQt ) )
    {
        return 1;
    }

    QImageWriter writer;
    writer.setFormat("png");
    writer.setFileName( outfile );
    if( !writer.write( *imageQt ) )
    {
        return 1;
    }

    return 0;
}

```

## 29.21 CreateARGBImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre

```

```

All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * http://www.w3.org/Graphics/PNG/inline-alpha.html
 * alphatest.png: PNG image data, 380 x 287, 8-bit/color RGBA, non-interlaced
 *
 * $ convert alphatest.png alphatest.rgba
 */

#include "gdcmImageReader.h"
#include "gdcmSequenceOfFragments.h"
#include "gdcmSystem.h"
#include "gdcmImageWriter.h"

#include <iostream>
#include <fstream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.rgb output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    size_t len = gdcm::System::FileSize(filename);
    std::ifstream is(filename, std::ios::binary);

    char * buf = new char[len];
    is.read(buf, len);

    gdcm::ImageWriter writer;
    gdcm::Image &image = writer.GetImage();
    image.SetNumberOfDimensions( 2 );
    unsigned int dims[3] = {};
    dims[0] = 380;
    dims[1] = 287;
    image.SetDimensions( dims );
    gdcm::PixelFormat pf = gdcm::PixelFormat::UINT8;
    pf.SetSamplesPerPixel( 4 );
    image.SetPixelFormat( pf );
    gdcm::PhotometricInterpretation pi =
        gdcm::PhotometricInterpretation::ARGB;
    image.SetPhotometricInterpretation( pi );
    image.SetTransferSyntax(
        gdcm::TransferSyntax::ExplicitVRLittleEndian );

    gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
    pixeldata.SetByteValue( buf, (uint32_t)len );
    image.SetDataElement( pixeldata );

    writer.SetFileName( outfile );
    if( !writer.Write() )
    {
        return 1;
    }
    delete[] buf;

    return 0;
}

```

## 29.22 CreateCMYKImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre

```

All rights reserved.  
See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even  
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR  
PURPOSE. See the above copyright notice for more information.

```
=====*/
/*
 * http://www.w3.org/Graphics/PNG/inline-alpha.html
 * alphatest.png: PNG image data, 380 x 287, 8-bit/color RGBA, non-interlaced
 *
 * $ convert alphatest.png alphatest.cmyk
 */

#include "gdcmImageReader.h"
#include "gdcmSequenceOfFragments.h"
#include "gdcmSystem.h"
#include "gdcmImageWriter.h"

#include <iostream>
#include <fstream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.cmyk output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    size_t len = gdcm::System::FileSize(filename);
    std::ifstream is(filename, std::ios::binary);

    char * buf = new char[len];
    is.read(buf, len);

    gdcm::ImageWriter writer;
    gdcm::Image &image = writer.GetImage();
    image.SetNumberOfDimensions( 2 );
    unsigned int dims[3] = {};
    dims[0] = 380;
    dims[1] = 287;
    image.SetDimensions( dims );
    gdcm::PixelFormat pf = gdcm::PixelFormat::UINT8;
    pf.SetSamplesPerPixel( 4 );
    image.SetPixelFormat( pf );
    gdcm::PhotometricInterpretation pi =
        gdcm::PhotometricInterpretation::CMYK;
    image.SetPhotometricInterpretation( pi );
    image.SetTransferSyntax(
        gdcm::TransferSyntax::ExplicitVRLittleEndian );

    gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
    pixeldata.SetByteValue( buf, (uint32_t)len );
    image.SetDataElement( pixeldata );

    writer.SetFileName( outfile );
    if( !writer.Write() )
    {
        return 1;
    }
    delete[] buf;

    return 0;
}
```

## 29.23 CreateFakeRTDOSE.cxx

```
/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
```

See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the above copyright notice for more information.

```

=====*/
#include "vtkGDCMImageWriter.h"
#include "vtkImageReader.h"
#include "vtkImageCast.h"
#include "vtkImageData.h"
#include "vtkPointData.h"
#include "vtkDataArray.h"
#include "vtkMedicalImageProperties.h"

#include "gdcmTrace.h"
#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmAttribute.h"

/*
 * Minimal example to create a fake RTDOSE file. The data contains a sphere
 * just for testing.
 * The vtkMedicalImageProperties is not properly filled, but only contains a
 * single field which is required to set the proper SOP Class
 */
int main(int, char *[])
{
    //gdcm::Trace::DebugOn();

    const vtkIdType xSize = 512;
    const vtkIdType ySize = 512;
    const vtkIdType zSize = 512;

    vtkImageData *image = vtkImageData::New();
    image->SetDimensions(xSize,ySize,zSize);
    image->SetOrigin(-350.684,350.0,890.76);
    image->SetSpacing(5.4688,-5.4688,-3.27);
    #if VTK_MAJOR_VERSION <= 5
        image->SetNumberOfScalarComponents(1);
        image->SetScalarTypeToDouble();
    #else
        image->AllocateScalars(VTK_DOUBLE ,1);
    #endif

    double pt[3];
    for( int z = 0; z < zSize; ++z )
        for( int y = 0; y < ySize; ++y )
            for( int x = 0; x < xSize; ++x )
            {
                pt[0] = x;
                pt[1] = y;
                pt[2] = z;
                pt[0] -= xSize / 2;
                pt[1] -= ySize / 2;
                pt[2] -= zSize / 2;
                pt[0] /= xSize / 2;
                pt[1] /= ySize / 2;
                pt[2] /= zSize / 2;
                const double unit = pt[0] * pt[0] + pt[1] * pt[1] + pt[2] * pt[2];
                const double inval = unit <= 1. ? (3 * unit + 7) : 0.; // just for fun => max == 10.
                double* pixel= static_cast<double*>(image->GetScalarPointer(x,y,z));
                pixel[0] = inval;
            }

    vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
    writer->SetFileDimensionality( 3 );
    writer->SetFileName( "rtdose.dcm" );
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputData( image );
    #else
        writer->SetInput( image );
    #endif
    writer->GetMedicalImageProperties()->SetSliceThickness("1.5");
    writer->GetMedicalImageProperties()->AddUserDefinedValue( "Dose Units", "GY");
    writer->GetMedicalImageProperties()->AddUserDefinedValue( "Dose Summation Type", "PLAN");
    writer->GetMedicalImageProperties()->AddUserDefinedValue( "Dose Type", "PHYSICAL");
    writer->GetMedicalImageProperties()->AddUserDefinedValue( "Frame of Reference UID", "
        1.3.12.2.1107.5.6.1.68100.30270111041215391275000000001");
    writer->GetMedicalImageProperties()->SetModality( "RTDOSE" );

```

```

writer->SetScale( 0.0042 ); // why not
writer->Write();

image->Delete();
writer->Delete();

// BEGIN HACK
// In GDCM version 2.4.3 and before, the following tag was missing which caused issue with some RTDose
// software:

// Open the DICOM file that was temporarily created. This will allows me to used
// GDCM to append specific tags that allows the RTDOSE to be associated with the
// relevant CT images.
gdcmm::Reader reader2;
reader2.SetFileName("rtdose.dcm" );
reader2.Read();
gdcmm::File &file = reader2.GetFile();
gdcmm::DataSet &ds = file.GetDataSet();

// Required by some software and not automagically added by GDCM in old version
gdcmm::Attribute<0x0028,0x0009> framePointer;
framePointer.SetNumberOfValues(1);
framePointer.SetValue( gdcmm::Tag(0x3004,0x000C) );
ds.Replace( framePointer.GetAsDataElement() );

gdcmm::Writer writer2;
writer2.CheckFileMetaInformationOff();
writer2.SetFileName("rtdose2.dcm");
writer2.SetFile( file );
writer2.Write();
// END HACK

return 0;
}

```

## 29.24 CreateJPIPDataSet.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example was created during the GSOC 2011 project for
 * JPIP
 */
#include "gdcmmAnonymizer.h"
#include "gdcmmWriter.h"
#include "gdcmmUIDGenerator.h"
#include "gdcmmFile.h"
#include "gdcmmTag.h"
#include "gdcmmSystem.h"
#include "gdcmmAttribute.h"

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " output.dcm" << std::endl;
        return 1;
    }
    const char *outfilename = argv[1];

    gdcmm::Writer w;
    gdcmm::File &file = w.GetFile();
    gdcmm::DataSet &ds = file.GetDataSet();
    //w.SetCheckFileMetaInformation( true );
    w.SetFileName( outfile );
}

```



```

file.GetHeader().SetDataSetTransferSyntax(
    gdcm::TransferSyntax::JPIPReferenced );

gdcm::Anonymizer anon;
anon.SetFile( file );

gdcm::MediaStorage ms =
    gdcm::MediaStorage::SecondaryCaptureImageStorage;

gdcm::UIDGenerator gen;
anon.Replace( gdcm::Tag(0x0008,0x16), ms.GetString() );
std::cout << ms.GetString() << std::endl;
anon.Replace( gdcm::Tag(0x0008,0x18), gen.Generate() );
//
anon.Replace( gdcm::Tag(0x0010,0x10), "JPIP^EXAMPLE" );
anon.Replace( gdcm::Tag(0x0010,0x20), "012345" );
anon.Empty( gdcm::Tag(0x0010,0x30) );
anon.Empty( gdcm::Tag(0x0010,0x40) );
anon.Empty( gdcm::Tag(0x0008,0x20) );
anon.Empty( gdcm::Tag(0x0008,0x30) );
anon.Empty( gdcm::Tag(0x0008,0x90) );
anon.Empty( gdcm::Tag(0x0020,0x10) );
anon.Empty( gdcm::Tag(0x0020,0x11) );
anon.Empty( gdcm::Tag(0x0008,0x50) );
anon.Empty( gdcm::Tag(0x0020,0x0013) );
anon.Replace( gdcm::Tag(0x0020,0xd), gen.Generate() );
anon.Replace( gdcm::Tag(0x0020,0xe), gen.Generate() );
anon.Replace( gdcm::Tag(0x0008,0x64), "WSD " );
anon.Replace( gdcm::Tag(0x0008,0x60), "OT" );

gdcm::Attribute<0x0028,0x7FE0> at;
at.SetValue( "http://dicom.example.com/jpipserver.cgi?target=img.jp2" );
ds.Insert( at.GetAsDataElement() );

// Need to retrieve the PixelFormat information from the given file

if (!w.Write() )
{
    std::cerr << "Could not write: " << outfilename << std::endl;
    return 1;
}

return 0;
}

```

## 29.25 CreateRAWStorage.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 ""
16 <uid value="1.2.840.10008.5.1.4.1.1.66" name="Raw Data Storage" type="SOP Class" part="PS 3.4" retired=
  "false"/>
17 ""
18
19 import gdcm
20 import sys,os
21
22 if __name__ == "__main__":
23     r = gdcm.Reader()
24     # Will require Testing...
25     dataroot = gdcm.Testing.GetDataRoot()
26     filename = os.path.join( dataroot, '012345.002.050.dcm' )
27     r.SetFileName( filename )
28     r.Read()

```

```
29  f = r.GetFile()
30  ds = f.GetDataSet()
31
32  uid = "1.2.840.10008.5.1.4.1.1.66"
33  # f = gdcM.File()
34  # ds = f.GetDataSet()
35  de = gdcM.DataElement( gdcM.Tag(0x0008,0x0016) )
36  de.SetByteValue( uid, gdcM.VL(len(uid)) )
37  vr = gdcM.VR( gdcM.VR.UI )
38  de.SetVR( vr )
39  ds.Replace( de )
40
41  ano = gdcM.Anonymizer()
42  ano.SetFile( r.GetFile() )
43  ano.RemovePrivateTags()
44  ano.RemoveGroupLength()
45  taglist = [
46  gdcM.Tag(0x0008,0x0008),
47  gdcM.Tag(0x0008,0x0022),
48  gdcM.Tag(0x0008,0x0032),
49  gdcM.Tag(0x0008,0x2111),
50  gdcM.Tag(0x0008,0x1150),
51  gdcM.Tag(0x0008,0x1155),
52  gdcM.Tag(0x0008,0x0100),
53  gdcM.Tag(0x0008,0x0102),
54  gdcM.Tag(0x0008,0x0104),
55  gdcM.Tag(0x0040,0xa170),
56  gdcM.Tag(0x0008,0x2112),
57  gdcM.Tag(0x0008,0x0100),
58  gdcM.Tag(0x0008,0x0102),
59  gdcM.Tag(0x0008,0x0104),
60  gdcM.Tag(0x0008,0x9215),
61  gdcM.Tag(0x0018,0x0010),
62  gdcM.Tag(0x0018,0x0022),
63  gdcM.Tag(0x0018,0x0050),
64  gdcM.Tag(0x0018,0x0060),
65  gdcM.Tag(0x0018,0x0088),
66  gdcM.Tag(0x0018,0x0090),
67  gdcM.Tag(0x0018,0x1040),
68  gdcM.Tag(0x0018,0x1100),
69  gdcM.Tag(0x0018,0x1110),
70  gdcM.Tag(0x0018,0x1111),
71  gdcM.Tag(0x0018,0x1120),
72  gdcM.Tag(0x0018,0x1130),
73  gdcM.Tag(0x0018,0x1150),
74  gdcM.Tag(0x0018,0x1151),
75  gdcM.Tag(0x0018,0x1152),
76  gdcM.Tag(0x0018,0x1160),
77  gdcM.Tag(0x0018,0x1190),
78  gdcM.Tag(0x0018,0x1210),
79  gdcM.Tag(0x0020,0x0012),
80  gdcM.Tag(0x0020,0x0032),
81  gdcM.Tag(0x0020,0x0037),
82  gdcM.Tag(0x0020,0x1041),
83  gdcM.Tag(0x0020,0x4000),
84  gdcM.Tag(0x0028,0x0002),
85  gdcM.Tag(0x0028,0x0004),
86  gdcM.Tag(0x0028,0x0010),
87  gdcM.Tag(0x0028,0x0011),
88  gdcM.Tag(0x0028,0x0030),
89  gdcM.Tag(0x0028,0x0100),
90  gdcM.Tag(0x0028,0x0101),
91  gdcM.Tag(0x0028,0x0102),
92  gdcM.Tag(0x0028,0x0103),
93  gdcM.Tag(0x0028,0x1052),
94  gdcM.Tag(0x0028,0x1053),
95  gdcM.Tag(0x0028,0x2110),
96  gdcM.Tag(0x0028,0x2112),
97  gdcM.Tag(0x7fe0,0x0010),
98  gdcM.Tag(0x0018,0x0020),
99  gdcM.Tag(0x0018,0x0021),
100 gdcM.Tag(0x0018,0x0023),
101 gdcM.Tag(0x0018,0x0025),
102 gdcM.Tag(0x0018,0x0080),
103 gdcM.Tag(0x0018,0x0081),
104 gdcM.Tag(0x0018,0x0083),
105 gdcM.Tag(0x0018,0x0084),
106 gdcM.Tag(0x0018,0x0085),
107 gdcM.Tag(0x0018,0x0086),
108 gdcM.Tag(0x0018,0x0087),
109 gdcM.Tag(0x0018,0x0091),
```

```

110     gdcM.Tag(0x0018,0x0093),
111     gdcM.Tag(0x0018,0x0094),
112     gdcM.Tag(0x0018,0x0095),
113     gdcM.Tag(0x0018,0x1088),
114     gdcM.Tag(0x0018,0x1090),
115     gdcM.Tag(0x0018,0x1094),
116     gdcM.Tag(0x0018,0x1250),
117     gdcM.Tag(0x0018,0x1251),
118     gdcM.Tag(0x0018,0x1310),
119     gdcM.Tag(0x0018,0x1312),
120     gdcM.Tag(0x0018,0x1314),
121     gdcM.Tag(0x0018,0x1315),
122     gdcM.Tag(0x0018,0x1316),
123     gdcM.Tag(0x0020,0x0110),
124     gdcM.Tag(0x0028,0x0120),
125     gdcM.Tag(0x0028,0x1050),
126     gdcM.Tag(0x0028,0x1051)
127 ]
128 for tag in taglist:
129     #print tag
130     ano.Remove( tag )
131
132     # special handling
133     gen = gdcM.UIDGenerator()
134     ano.Replace( gdcM.Tag(0x0008,0x9123), gen.Generate() )
135     #ano.Empty( gdcM.Tag(0x0040,0x0555) )
136
137
138 #
139 # uid = gen.Generate()
140 # de.SetTag( gdcM.Tag(0x0008,0x0018) )
141 # de.SetByteValue( uid, gdcM.VL(len(uid)) )
142 # ds.Insert( de )
143
144 # init FMI now:
145 #fmi = f.GetHeader()
146 #ts = gdcM.TransferSyntax()
147 #print ts
148 #fmi.SetDataSetTransferSyntax( ts ) # default
149 #print fmi.GetDataSetTransferSyntax()
150 #de.SetTag( gdcM.Tag(0x0002,0x0010) )
151 #uid = "1.2.840.10008.1.2"
152 #de.SetByteValue( uid, gdcM.VL(len(uid)) )
153 #fmi.Insert( de )
154 # f.SetHeader( r.GetFile().GetHeader() )
155
156 writer = gdcM.Writer()
157 writer.SetFile( ano.GetFile() )
158 writer.SetFileName( "rawstorage.dcm" );
159 writer.Write()

```

## 29.26 csa2img.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * I do not know what the format is, just guessing from info found on the net:
 *
 * http://atonal.ucdavis.edu/matlab/fmri/spm5/spm_dicom_convert.m
 *
 * This example is an attempt at understanding the format used by SIEMENS
 * their "SIEMENS CSA NON-IMAGE" DICOM file (1.3.12.2.1107.5.9.1)
 *
 * Everything done in this code is for the sole purpose of writing interoperable
 * software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
 * If you believe anything in this code violates any law or any of your rights,

```

```

* please contact us (gdcm-developers@lists.sourceforge.net) so that we can
* find a solution.
*
*/
#include "gdcmReader.h"
#include "gdcmImageReader.h"
#include "gdcmImageWriter.h"
#include "gdcmCSAHeader.h"
#include "gdcmAttribute.h"
#include "gdcmPrivateTag.h"

#include <math.h>

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    // gdcmDataExtra/gdcmNonImageData/exCSA_Non-Image_Storage.dcm
    // PHANTOM.MR.CARDIO_COEUR_S_QUENCE_DE_REP_RANGE.9.257.2008.03.20.14.53.25.578125.43151705.IMA
    const char *filename = argv[1];

    gdcm::Reader reader; // Do not use ImageReader
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }

    gdcm::CSAHeader csa;
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    const gdcm::PrivateTag &t1 = csa.GetCSAImageHeaderInfoTag();
    //std::cout << t1 << std::endl;
    //const gdcm::PrivateTag &t2 = csa.GetCSASeriesHeaderInfoTag();

    if( ds.FindDataElement( t1 ) )
    {
        csa.LoadFromDataElement( ds.GetDataElement( t1 ) );
        csa.Print( std::cout );
    }

    int dims[2] = {};
    if( csa.FindCSAElementByName( "Columns" ) )
    {
        const gdcm::CSAElement &crael = csa.GetCSAElementByName( "Columns" )
        ;
        std::cout << crael << std::endl;
        //const gdcm::ByteValue *bv = crael.GetByteValue();
        gdcm::Element<gdcm::VR::IS, gdcm::VM::VM1> el;
        el.Set( crael.GetValue() );
        dims[0] = el.GetValue();
        std::cout << "Columns:" << el.GetValue() << std::endl;
    }

    if( csa.FindCSAElementByName( "Rows" ) )
    {
        const gdcm::CSAElement &crael2 = csa.GetCSAElementByName( "Rows" );
        std::cout << crael2 << std::endl;
        gdcm::Element<gdcm::VR::IS, gdcm::VM::VM1> el2;
        el2.Set( crael2.GetValue() );
        dims[1] = el2.GetValue();
        std::cout << "Rows:" << el2.GetValue() << std::endl;
    }

    double spacing[2] = { 1. , 1. };
    bool spacingfound = false;
    if( csa.FindCSAElementByName( "PixelSpacing" ) )
    {
        const gdcm::CSAElement &crael3 = csa.GetCSAElementByName( "
        PixelSpacing" );
        if( !crael3.IsEmpty() )
        {
            std::cout << crael3 << std::endl;
            gdcm::Element<gdcm::VR::DS, gdcm::VM::VM2> el3;
            el3.Set( crael3.GetValue() );
            spacing[0] = el3.GetValue(0);
            spacing[1] = el3.GetValue(1);
            std::cout << "PixelSpacing:" << el3.GetValue() << "," << el3.
            GetValue(1) << std::endl;
            spacingfound = true;
        }
    }
}

```

```

if( !spacingfound )
{
    std::cerr << "Problem with PixelSpacing" << std::endl;
    //return 1;
}
if( !dims[0] || !dims[1] )
{
    std::cerr << "Problem with dims" << std::endl;
    return 1;
}

gdcm::ImageWriter writer;

gdcm::Image &image = writer.GetImage();
image.SetNumberOfDimensions( 2 ); // good default
image.SetDimension(0, dims[0] );
image.SetDimension(1, dims[1] );
image.SetSpacing(0, spacing[0] );
image.SetSpacing(1, spacing[1] );
gdcm::PixelFormat pixeltype = gdcm::PixelFormat::INT16; //
    bytewidth = spm_type('int16','bits')/8;

//unsigned long l = image.GetBufferLength();
//const int p = l / (dims[0] * dims[1]);

//image.SetNumberOfDimensions( 3 );
//image.SetDimension(2, p / pixeltype.GetPixelSize() );

gdcm::PhotometricInterpretation pi;
pi = gdcm::PhotometricInterpretation::MONOCHROME2;
//pixeltype.SetSamplesPerPixel( );
image.SetPhotometricInterpretation( pi );
image.SetPixelFormat( pixeltype );
//image.SetIntercept( inputimage.GetIntercept() );
//image.SetSlope( inputimage.GetSlope() );

//gdcm::DataElement pixeldata( gdcm::Tag(0x7fe1,0x1010) );
//pixeldata.SetByteValue( &outbuf[0], outbuf.size() );
gdcm::PrivateTag csanonimaget(0x7fe1,0x10,"SIEMENS CSA NON-IMAGE");
const gdcm::DataElement &pixeldata = ds.GetDataElement( csanonimaget );
image.SetDataElement( pixeldata );

std::string outfilename = "outcsa.dcm";
//writer.SetFile( reader.GetFile() );
writer.SetFileName( outfilename.c_str() );
if( !writer.Write() )
{
    std::cerr << "could not write: " << outfilename << std::endl;
    return 1;
}

return 0;
}

```

## 29.27 CStoreQtProgress.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This small example show how one can use the virtual function
 * mechanism of the SimpleSubjectWatcher class to redirect progress
 * report to a custom Qt classes
 *
 * http://doc.qt.nokia.com/latest/qprogressdialog.html

```

```

*
* Usage:
* CStoreQtProgress dicom.example.com 11112 gdcmlData/MR_Spectroscopy_SIEMENS_OF.dcm
*
*/

#include "gdcmlServiceClassUser.h"
#include "gdcmlSimpleSubjectWatcher.h"
#include "gdcmlProgressEvent.h"
#include "gdcmlDirectory.h"
#include "gdcmlPresentationContextGenerator.h"

#include <QApplication>
#include <QProgressDialog>
#include <QVBoxLayout>

namespace gdcml {
/*
 * This class is a little more complicated than what this example demonstrate
 * This watcher is capable of handling nested progress. Since the Progress
 * grows from [0 to 1] on a per file basis and we only have one instance of a
 * watcher per association, we need some calculation to compute the global
 * (total) progress
 * In fact we simply divide the per-file progress by the number of files.
 *
 * This QtWatcher class will then update the progress bar according to the
 * progress.
 */
class MyQtWatcher : public SimpleSubjectWatcher
{
    size_t nfiles;
    double progress;
    size_t index;
    double refprogress;
    QWidget* win;
    QProgressDialog* qtprogress;
public:
    MyQtWatcher(Subject * s, const char *comment = "", QWidget *w = NULL, QProgressDialog* p = NULL, size_t n
        = 1):
        SimpleSubjectWatcher(s,comment),nfiles(n),progress(0),index(0),refprogress(0),win(w),qtprogress(p) {}
    void ShowIteration()
    {
        index++;
        assert( index <= nfiles );
        // update refprogress (we are moving to the next file)
        refprogress = progress;
    }
    void ShowProgress(Subject *, const Event &evt)
    {
        // Retrieve the ProgressEvent:
        const ProgressEvent &pe = dynamic_cast<const ProgressEvent&>(evt);
        // compute global progress:
        progress = refprogress + (1. / (double)nfiles ) * pe.GetProgress();
        // Print Global and local progress to stdout:
        std::cout << "Global Progress: " << progress << " per file progress " << pe.GetProgress() << std::endl;
        //set progress value in the QtProgress bar
        int i = (int)(progress * 100 + 0.5); // round to next int
        qtprogress->setValue(i);
        win->show();
    }
    virtual void ShowDataSet(Subject *caller, const Event &evt)
    {
        (void)caller;
        (void)evt;
    }
};
} // end namespace gdcml

int main(int argc, char *argv[])
{
    if( argc < 4 )
    {
        std::cerr << argv[0] << " remote_server port filename" << std::endl;
        return 1;
    }
    QApplication a(argc, argv);

    std::ostringstream error_log;
    gdcml::Trace::SetErrorStream( error_log );

    const char *remote = argv[1];

```

```

int portno = atoi(argv[2]);
const char *filename = argv[3];

QVBoxLayout* layout = new QVBoxLayout;
QWidget* win = new QWidget;

QProgressDialog* progress = new QProgressDialog("Sending data...", "Cancel", 0, 100);
progress->setWindowModality(Qt::WindowModal);

layout->addWidget(progress, Qt::AlignCenter);
win->setLayout(layout);

gdcm::SmartPointer<gdcm::ServiceClassUser> scup = new
    gdcm::ServiceClassUser;
gdcm::ServiceClassUser &scu = *scup;
//gdcm::SimpleSubjectWatcher w( &scu, "TestServiceClassUser" );
// let's use a more complicated progress reported in this example
gdcm::MyQtWatcher w( &scu, "QtWatcher", win, progress );

scu.SetHostname( remote );
scu.SetPort( (uint16_t)portno );
scu.SetTimeout( 1000 );
scu.SetCalledAETitle( "GDCM_STORE" );

if( !scu.InitializeConnection() )
{
    std::cerr << "Could not InitializeConnection" << std::endl;
    return 1;
}

gdcm::Directory::FileNamesType filenames;
filenames.push_back( filename );

// setup the PC(s) based on the filenames:
gdcm::PresentationContextGenerator generator;
if( !generator.GenerateFromFileNames(filenames) )
{
    std::cerr << "Could not GenerateFromFileNames" << std::endl;
    return 1;
}

// Setup PresentationContext(s)
scu.SetPresentationContexts( generator.
    GetPresentationContexts() );

// Start ASSOCIATION
if( !scu.StartAssociation() )
{
    std::cerr << "Could not Start" << std::endl;
    return 1;
}

// Send C-STORE
if( !scu.SendStore( filename ) )
{
    std::cerr << "Could not Store" << std::endl;
    std::cerr << "Error log is:" << std::endl;
    std::cerr << error_log.str() << std::endl;
    return 1;
}

// Stop ASSOCIATION
if( !scu.StopAssociation() )
{
    std::cerr << "Could not Stop" << std::endl;
    return 1;
}

win->show();

return a.exec();
}

```

## 29.28 DecompressImage.cs

This is a C# example on how to use Image

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/DecompressImage.exe gdcmData/012345.002.050.dcm decompress.dcm
 */
using System;
using gdcm;

public class DecompressImage
{
    public static int Main(string[] args)
    {
        string file1 = args[0];
        string file2 = args[1];
        ImageReader reader = new ImageReader();
        reader.SetFileName( file1 );
        bool ret = reader.Read();
        if( !ret )
        {
            return 1;
        }

        Image image = new Image();
        Image ir = reader.GetImage();

        image.SetNumberOfDimensions( ir.GetNumberOfDimensions() );

        //Just for fun:
        //int dircos = ir.GetDirectionCosines();
        //t = gdcm.Orientation.GetType(dircos);
        //int l = gdcm.Orientation.GetLabel(t);
        //System.Console.WriteLine( "Orientation label:" + l );

        // Set the dimensions,
        // 1. either one at a time
        //image.SetDimension(0, ir.GetDimension(0) );
        //image.SetDimension(1, ir.GetDimension(1) );

        // 2. the array at once
        uint[] dims = {0, 0};
        // Just for fun let's invert the dimensions:
        dims[0] = ir.GetDimension(1);
        dims[1] = ir.GetDimension(0);
        ir.SetDimensions( dims );

        PixelFormat pixeltype = ir.GetPixelFormat();
        image.SetPixelFormat( pixeltype );

        PhotometricInterpretation pi = ir.GetPhotometricInterpretation();
        image.SetPhotometricInterpretation( pi );

        DataElement pixeldata = new DataElement( new Tag(0x7fe0,0x0010) );
        byte[] str1 = new byte[ ir.GetBufferLength()];
        ir.GetBuffer( str1 );
        //System.Console.WriteLine( ir.GetBufferLength() );
        pixeldata.SetByteValue( str1, new VL( (uint)str1.Length ) );
        //image.SetDataElement( pixeldata );
        ir.SetDataElement( pixeldata );

        ImageWriter writer = new ImageWriter();
        writer.SetFileName( file2 );
        writer.SetFile( reader.GetFile() );
        writer.SetImage( ir );
        ret = writer.Write();
        if( !ret )
        {

```



```

        return 1;
    }

    return 0;
}

```

## 29.29 DecompressImage.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * This example will take in a DICOM file, and tries to decompress it (actually write it
 * as ImplicitVRLittleEndian Transfer Syntax).
 *
 * Compilation:
 * $ CLASSPATH=gdcm.jar javac ../../gdcm/Examples/Java/DecompressImage.java -d .
 *
 * Usage:
 * $ LD_LIBRARY_PATH=. CLASSPATH=gdcm.jar:. java DecompressImage gdcmData/012345.002.050.dcm out.dcm
 */
import gdcm.*;

public class DecompressImage
{
    public static void main(String[] args) throws Exception
    {
        String file1 = args[0];
        String file2 = args[1];
        ImageReader reader = new ImageReader();
        reader.SetFileName( file1 );
        boolean ret = reader.Read();
        if( !ret )
        {
            throw new Exception("Could not read: " + file1 );
        }

        ImageChangeTransferSyntax change = new ImageChangeTransferSyntax();
        change.SetTransferSyntax( new TransferSyntax(TransferSyntax.TSType.ImplicitVRLittleEndian) );
        change.SetInput( reader.GetImage() );
        if( !change.Change() )
        {
            throw new Exception("Could not change: " + file1 );
        }

        Image out = change.GetOutput();
        System.out.println( out.toString() );

        // Set the Source Application Entity Title
        FileMetaInformation.SetSourceApplicationEntityTitle( "Just For Fun" );

        ImageWriter writer = new ImageWriter();
        writer.SetFileName( file2 );
        writer.SetFile( reader.GetFile() );
        writer.SetImage( out );
        ret = writer.Write();
        if( !ret )
        {
            throw new Exception("Could not write: " + file2 );
        }
    }
}

```

## 29.30 DecompressImage.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #       This software is distributed WITHOUT ANY WARRANTY; without even
10 #       the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #       PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18   python DecompressImage.py gdcmData/012345.002.050.dcm decompress.dcm
19 """
20
21 import gdcm
22 import sys
23
24 if __name__ == "__main__":
25
26     file1 = sys.argv[1]
27     file2 = sys.argv[2]
28
29     r = gdcm.ImageReader()
30     r.SetFileName( file1 )
31     if not r.Read():
32         sys.exit(1)
33
34     image = gdcm.Image()
35     ir = r.GetImage()
36
37     image.SetNumberOfDimensions( ir.GetNumberOfDimensions() );
38     dims = ir.GetDimensions();
39     print ir.GetDimension(0);
40     print ir.GetDimension(1);
41     print "Dims:", dims
42
43     # Just for fun:
44     dircos = ir.GetDirectionCosines()
45     t = gdcm.Orientation.GetType(dircos)
46     l = gdcm.Orientation.GetLabel(t)
47     print "Orientation label:", l
48
49     image.SetDimension(0, ir.GetDimension(0) );
50     image.SetDimension(1, ir.GetDimension(1) );
51
52     pixeltype = ir.GetPixelFormat();
53     image.SetPixelFormat( pixeltype );
54
55     pi = ir.GetPhotometricInterpretation();
56     image.SetPhotometricInterpretation( pi );
57
58     pixeldata = gdcm.DataElement( gdcm.Tag(0x7fe0,0x0010) )
59     str1 = ir.GetBuffer()
60     #print ir.GetBufferLength()
61     pixeldata.SetByteValue( str1, gdcm.VL( len(str1) ) )
62     image.SetDataElement( pixeldata )
63
64     w = gdcm.ImageWriter()
65     w.SetFileName( file2 )
66     w.SetFile( r.GetFile() )
67     w.SetImage( image )
68     if not w.Write():
69         sys.exit(1)

```

## 29.31 DecompressImageMultiframe.cs

```

/*=====

```

```

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
$ gdcminfo ~/Desktop/angiogram-06.dcm
MediaStorage is 1.2.840.10008.5.1.4.1.1.12.1 [X-Ray Angiographic Image Storage]
TransferSyntax is 1.2.840.10008.1.2.4.50 [JPEG Baseline (Process 1): Default Transfer Syntax for Lossy JPEG
      8 Bit Image Compression]
NumberOfDimensions: 3
Dimensions: (512,512,355)
Origin: (0,0,0)
Spacing: (1,1,40)
DirectionCosines: (1,0,0,0,1,0)
Rescale Intercept/Slope: (0,1)
SamplesPerPixel      :1
BitsAllocated        :8
BitsStored           :8
HighBit              :7
PixelRepresentation:0
ScalarType found     :UINT8
PhotometricInterpretation: MONOCHROME2
PlanarConfiguration: 0
TransferSyntax: 1.2.840.10008.1.2.4.50
Orientation Label: AXIAL
*/

/*
* Description:
*
* Assume we have a file angiogram-06.dcm as described above.
* the following program will decompress directly from the extracted jpeg stream.
*
* First step extract the jpeg stream (but not the Basic Offset Table):
*
* $ gdcmrw -i angiogram-06.dcm -o /tmp/output/chris --split-frags --pattern %d.jpg
*
* Check that indeed there are 355 files, while there are 356 fragments in the original DICOM file, since
* gdcmrw always skip the first fragment (Basic Offset Table).
*
* Now from those individual jpeg stream, recreate a fake gdcm.DataElement...
*
* Usage:
*
* $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
* $ mono ./bin/DecompressImageMultiframe.exe /tmp/output
*/
using System;
using gdcm;

public class DecompressImageMultiframe
{
    public static int Main(string[] args)
    {
        string directory = args[0];
        gdcm.Directory dir = new gdcm.Directory();
        uint nfiles = dir.Load(directory);
        //System.Console.WriteLine(dir.toString());
        gdcm.FilenamesType filenames = dir.GetFilesNames();

        Image image = new Image();
        image.SetNumberOfDimensions( 3 ); // important for now
        DataElement pixeldata = new DataElement( new gdcm.Tag(0x7fe0,0x0010) );

        // Create a new SequenceOfFragments C++ object, store it as a SmartPointer :
        SmartPtrFrag sq = SequenceOfFragments.New();

        // Yeah, the file are not guarantee to be in order, please adapt...
        for(uint i = 0; i < nfiles; ++i)
        {
            System.Console.WriteLine( filenames[(int)i] );
            string file = filenames[(int)i];
            System.IO.FileStream infile =

```

```

        new System.IO.FileStream(file, System.IO.FileMode.Open, System.IO.FileAccess.Read);
        uint fsize = gdcm.PosixEmulation.FileSize(file);

        byte[] jstream = new byte[fsize];
        infile.Read(jstream, 0 , jstream.Length);

        Fragment frag = new Fragment();
        frag.SetByteValue( jstream, new gdcm.VL( (uint)jstream.Length) );
        sq.AddFragment( frag );
    }

    // Pass by reference:
    pixeldata.SetValue( sq.__ref__() );

    // insert:
    image.SetDataElement( pixeldata );

    // JPEG use YBR to achieve better compression ratio by default (not RGB)
    // FIXME hardcoded:
    PhotometricInterpretation pi = new PhotometricInterpretation( PhotometricInterpretation.PIType.
        MONOCHROME2 );
    image.SetPhotometricInterpretation( pi );
    // FIXME hardcoded:
    PixelFormat pixeltype = new PixelFormat(1,8,8,7);
    image.SetPixelFormat( pixeltype );

    // FIXME hardcoded:
    image.SetTransferSyntax( new TransferSyntax( TransferSyntax.TSType.JPEGLosslessProcess14_1 ) );
    image.SetDimension(0, 512);
    image.SetDimension(1, 512);
    image.SetDimension(2, 355);

    // Decompress !
    byte[] decompressedData = new byte[(int)image.GetBufferLength()];
    image.GetBuffer(decompressedData);

    // Write out the decompressed bytes
    System.Console.WriteLine(image.toString());
    using (System.IO.Stream stream =
        System.IO.File.Open(@"tmp/dd.raw",
            System.IO.FileMode.Create))
    {
        System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
        writer.Write(decompressedData);
    }

    return 0;
}

```

## 29.32 DecompressJPEGFile.cs

This is a C# example on how to use [gdcm::SequenceOfFragments](#)

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/DecompressJPEGFile.exe somejpegfile.jpg
 */
using System;
using gdcm;

```

```

public class DecompressJPEGFile
{
    public static int Main(string[] args)
    {
        string file1 = args[0];
        System.IO.FileStream infile =
            new System.IO.FileStream(file1, System.IO.FileMode.Open, System.IO.FileAccess.Read);
        uint fsize = gdcm.PosixEmulation.FileSize(file1);

        byte[] jstream = new byte[fsize];
        infile.Read(jstream, 0, jstream.Length);

        Trace.DebugOn();
        Image image = new Image();
        image.SetNumberOfDimensions( 2 ); // important for now
        DataElement pixeldata = new DataElement( new gdcm.Tag(0x7fe0,0x0010) );

        // DO NOT set a ByteValue here, JPEG is a particular kind of encapsulated syntax
        // in which can one cannot use a simple byte array for storage. Instead, see
        // gdcm.SequenceOfFragments
        //pixeldata.SetByteValue( jstream, new gdcm.VL( (uint)jstream.Length ) );

        // Create a new SequenceOfFragments C++ object, store it as a SmartPointer :
        SmartPtrFrag sq = SequenceOfFragments.New();
        Fragment frag = new Fragment();
        frag.SetByteValue( jstream, new gdcm.VL( (uint)jstream.Length ) );
        // Single file => single fragment
        sq.AddFragment( frag );
        // Pass by reference:
        pixeldata.SetValue( sq.__ref__() );

        // insert:
        image.SetDataElement( pixeldata );

        // JPEG use YBR to achieve better compression ratio by default (not RGB)
        // FIXME hardcoded:
        PhotometricInterpretation pi = new PhotometricInterpretation( PhotometricInterpretation.PIType.YBR_FULL
        );
        image.SetPhotometricInterpretation( pi );
        // FIXME hardcoded:
        PixelFormat pixeltype = new PixelFormat(3,8,8,7);
        image.SetPixelFormat( pixeltype );

        // FIXME hardcoded:
        image.SetTransferSyntax( new TransferSyntax( TransferSyntax.TSType.JPEGLosslessProcess14_1 ) );
        image.SetDimension(0, 692);
        image.SetDimension(1, 721);

        // Decompress !
        byte[] decompressedData = new byte[(int)image.GetBufferLength()];
        image.GetBuffer(decompressedData);

        // Write out the decompressed bytes
        System.Console.WriteLine(image.toString());
        using (System.IO.Stream stream =
            System.IO.File.Open(@"tmp/dd.raw",
                System.IO.FileMode.Create))
        {
            System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
            writer.Write(decompressedData);
        }

        return 0;
    }
}

```

## 29.33 DecompressPixmap.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

```

This software is distributed WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the above copyright notice for more information.

```

=====*/

/*
 * This example will take in a DICOM file, and tries to decompress it (actually write it
 * as ImplicitVRLittleEndian Transfer Syntax).
 *
 * Compilation:
 * $ CLASSPATH=gdcmm.jar javac ../../gdcmm/Examples/Java/DecompressPixmap.java -d .
 *
 * Usage:
 * $ LD_LIBRARY_PATH=. CLASSPATH=gdcmm.jar:. java DecompressPixmap gdcmmData/012345.002.050.dcm out.dcm
 */
import gdcmm.*;

public class DecompressPixmap
{
    public static void main(String[] args) throws Exception
    {
        String file1 = args[0];
        String file2 = args[1];
        PixmapReader reader = new PixmapReader();
        reader.SetFileName( file1 );
        boolean ret = reader.Read();
        if( !ret )
        {
            throw new Exception("Could not read: " + file1 );
        }

        ImageChangeTransferSyntax change = new ImageChangeTransferSyntax();
        change.SetTransferSyntax( new TransferSyntax(TransferSyntax.TSType.ImplicitVRLittleEndian) );
        PixmapToPixmapFilter filter = (PixmapToPixmapFilter)change;
        filter.SetInput( reader.GetPixmap() );
        if( !change.Change() )
        {
            throw new Exception("Could not change: " + file1 );
        }

        // The following does not work in Java/swig 2.0.7
        // Pixmap p = ((PixmapToPixmapFilter)change).GetOutput();
        Pixmap p = change.GetOutputAsPixmap(); // be explicit
        //System.out.println( p.toString() );

        // Set the Source Application Entity Title
        FileMetaInformation.SetSourceApplicationEntityTitle( "Just For Fun" );

        PixmapWriter writer = new PixmapWriter();
        writer.SetFileName( file2 );
        writer.SetFile( reader.GetFile() );
        writer.SetImage( p );
        ret = writer.Write();
        if( !ret )
        {
            throw new Exception("Could not write: " + file2 );
        }
    }
}

```

## 29.34 DiffFile.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

```

```

=====*/
#include "gdcmReader.h"

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input1.dcm input2.dcm" << std::endl;
        return 1;
    }
    const char *filename1 = argv[1];
    const char *filename2 = argv[2];

    gdcm::Reader reader1;
    reader1.SetFileName( filename1 );
    if( !reader1.Read() )
    {
        return 1;
    }

    gdcm::Reader reader2;
    reader2.SetFileName( filename2 );
    if( !reader2.Read() )
    {
        return 1;
    }

    const gdcm::File &file1 = reader1.GetFile();
    const gdcm::File &file2 = reader2.GetFile();

    const gdcm::DataSet &ds1 = file1.GetDataSet();
    const gdcm::DataSet &ds2 = file2.GetDataSet();

    gdcm::DataSet::ConstIterator it1 = ds1.Begin();
    gdcm::DataSet::ConstIterator it2 = ds2.Begin();

    const gdcm::DataElement &de1 = *it1;
    const gdcm::DataElement &de2 = *it2;
    if( de1 == de2 )
    {
    }
    while( it1 != ds1.End() && it2 != ds2.End() && *it1 == *it2 )
    {
        ++it1;
        ++it2;
    }

    if( it1 != ds1.End() || it2 != ds2.End() )
    {
        std::cerr << "Problem with:" << std::endl;
        if( it1 != ds1.End() )
        {
            std::cerr << "ds1: " << *it1 << std::endl;
        }
        if( it2 != ds2.End() )
        {
            std::cerr << "ds2: " << *it2 << std::endl;
        }
        return 1;
    }

    return 0;
}

```

## 29.35 DiscriminateVolume.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

```

```

=====*/
#include "gdcmScanner.h"
#include "gdcmTesting.h"
#include "gdcmIPPSorter.h"
#include "gdcmDirectionCosines.h"
#include <cmath>

/*
 * The following example is a basic sorted which should work in generic cases.
 * It sort files based on:
 * Study Instance UID
 * Series Instance UID
 * Frame of Reference UID
 * Image Orientation (Patient)
 * Image Position (Patient) (Sorting based on IPP + IOP)
 */

namespace gdcm {
    const Tag t1(0x0020,0x000d); // Study Instance UID
    const Tag t2(0x0020,0x000e); // Series Instance UID
    const Tag t3(0x0020,0x0052); // Frame of Reference UID
    const Tag t4(0x0020,0x0037); // Image Orientation (Patient)

    class DiscriminateVolume
    {
    private:
        std::vector< Directory::FileNamesType > SortedFiles;
        std::vector< Directory::FileNamesType > UnsortedFiles;

        Directory::FileNamesType GetAllFileNamesFromTagToValue(
            Scanner const & s, Directory::FileNamesType const & filesubset, Tag const & t,
            const char *valueref)
        {
            Directory::FileNamesType theReturn;
            if( valueref )
            {
                size_t len = strlen( valueref );
                Directory::FileNamesType::const_iterator file = filesubset.begin();
                for(; file != filesubset.end(); ++file)
                {
                    const char *filename = file->c_str();
                    const char * value = s.GetValue(filename, t);
                    if( value && strncmp(value, valueref, len ) == 0 )
                    {
                        theReturn.push_back( filename );
                    }
                }
            }
            return theReturn;
        }
    };

    void ProcessAIOP(Scanner const & , Directory::FileNamesType const & subset, const
        char *iopval)
    {
        std::cout << "IOP: " << iopval << std::endl;
        IPPSorter ipp;
        ipp.SetComputeZSpacing( true );
        ipp.SetZSpacingTolerance( 1e-3 ); // ??
        bool b = ipp.Sort( subset );
        if( !b )
        {
            // If you reach here this means you need one more parameter to discriminiat this
            // series. Eg. T1 / T2 intertwined. Multiple Echo (0018,0081)
            std::cerr << "Failed to sort: " << subset.begin()->c_str() << std::endl;
            for(
                Directory::FileNamesType::const_iterator file = subset.begin();
                file != subset.end(); ++file)
            {
                std::cerr << *file << std::endl;
            }
            UnsortedFiles.push_back( subset );
            return ;
        }
        ipp.Print( std::cout );
        SortedFiles.push_back( ipp.GetFileNames() );
    }

    void ProcessAFrameOfRef(Scanner const & s, Directory::FileNamesType const & subset,
        const char * frameuid)
    {

```



```

// In this subset of files (belonging to same series), let's find those
// belonging to the same Frame ref UID:
Directory::FileNamesType files = GetAllFileNamesFromTagToValue(
    s, subset, t3, frameuid);

std::set< std::string > iopset;

for(
    Directory::FileNamesType::const_iterator file = files.begin();
    file != files.end(); ++file)
{
    //std::cout << *file << std::endl;
    const char * value = s.GetValue(file->c_str(), gdc::t4 );
    assert( value );
    iopset.insert( value );
}

size_t n = iopset.size();
if ( n == 0 )
{
    assert( files.empty() );
    return;
}

std::cout << "Frame of Ref: " << frameuid << std::endl;
if ( n == 1 )
{
    ProcessAIOP(s, files, iopset.begin()->c_str() );
}
else
{
    const char *f = files.begin()->c_str();
    std::cerr << "More than one IOP: " << f << std::endl;
    // Make sure that there is actually 'n' different IOP
    gdc::DirectionCosines ref;
    gdc::DirectionCosines dc;
    for(
        std::set< std::string >::const_iterator it = iopset.begin();
        it != iopset.end(); ++it )
    {
        ref.SetFromString( it->c_str() );
        for(
            Directory::FileNamesType::const_iterator file = files.begin();
            file != files.end(); ++file)
        {
            std::string value = s.GetValue(file->c_str(), gdc::t4 );
            if( value != it->c_str() )
            {
                dc.SetFromString( value.c_str() );
                const double crossdot = ref.CrossDot(dc);
                const double eps = std::fabs( 1. - crossdot );
                if( eps < 1e-6 )
                {
                    std::cerr << "Problem with IOP discrimination: " << file->c_str()
                        << " " << it->c_str() << std::endl;
                    return;
                }
            }
        }
    }
    // If we reach here this means there is actually 'n' different IOP
    for(
        std::set< std::string >::const_iterator it = iopset.begin();
        it != iopset.end(); ++it )
    {
        const char *iopvalue = it->c_str();
        Directory::FileNamesType iopfiles = GetAllFileNamesFromTagToValue(
            s, files, t4, iopvalue );
        ProcessAIOP(s, iopfiles, iopvalue );
    }
}

void ProcessASeries(Scanner const & s, const char * seriesuid)
{
    std::cout << "Series: " << seriesuid << std::endl;
    // let's find all files belonging to this series:
    Directory::FileNamesType seriesfiles = GetAllFileNamesFromTagToValue(
        s, s.GetFileNames(), t2, seriesuid);

    gdc::Scanner::ValuesType vt3 = s.GetValues(t3);
    for(

```

```

        gdcM::Scanner::ValuesType::const_iterator it = vt3.begin()
        ; it != vt3.end(); ++it )
    {
        ProcessAFrameOfRef(s, seriesfiles, it->c_str());
    }
}

void ProcessAStudy(Scanner const & s, const char * studyuid)
{
    std::cout << "Study: " << studyuid << std::endl;
    gdcM::Scanner::ValuesType vt2 = s.GetValues(t2);
    for(
        gdcM::Scanner::ValuesType::const_iterator it = vt2.begin()
        ; it != vt2.end(); ++it )
    {
        ProcessASeries(s, it->c_str());
    }
}
public:

void Print( std::ostream & os )
{
    os << "Sorted Files: " << std::endl;
    for(
        std::vector< Directory::FilenameType >::const_iterator it = SortedFiles.begin();
        it != SortedFiles.end(); ++it )
    {
        os << "Group: " << std::endl;
        for(
            Directory::FilenameType::const_iterator file = it->begin();
            file != it->end(); ++file)
        {
            os << *file << std::endl;
        }
    }
    os << "Unsorted Files: " << std::endl;
    for(
        std::vector< Directory::FilenameType >::const_iterator it = UnsortedFiles.begin();
        it != UnsortedFiles.end(); ++it )
    {
        os << "Group: " << std::endl;
        for(
            Directory::FilenameType::const_iterator file = it->begin();
            file != it->end(); ++file)
        {
            os << *file << std::endl;
        }
    }
}

std::vector< Directory::FilenameType > const & GetSortedFiles() const { return SortedFiles; }
std::vector< Directory::FilenameType > const & GetUnsortedFiles() const { return UnsortedFiles; }

void ProcessIntoVolume( Scanner const & s )
{
    gdcM::Scanner::ValuesType vt1 = s.GetValues( gdcM::t1 );
    for(
        gdcM::Scanner::ValuesType::const_iterator it = vt1.begin()
        ; it != vt1.end(); ++it )
    {
        ProcessAStudy( s, it->c_str() );
    }
}

};

} // namespace gdcM

int main(int argc, char *argv[])
{
    std::string dirl;
    if( argc < 2 )
    {
        const char *extradataroot = NULL;
#ifdef GDCM_BUILD_TESTING
        extradataroot = gdcM::Testing::GetDataExtraRoot();
#endif
        if( !extradataroot )
        {

```

```

        return 1;
    }
    dir1 = extradataroot;
    dir1 += "/gdcmSampleData/ForSeriesTesting/VariousIncidences/ST1";
}
else
{
    dir1 = argv[1];
}

gdcm::Directory d;
d.Load( dir1.c_str(), true ); // recursive !

gdcm::Scanner s;
s.AddTag( gdcm::t1 );
s.AddTag( gdcm::t2 );
s.AddTag( gdcm::t3 );
s.AddTag( gdcm::t4 );
bool b = s.Scan( d.GetFilesNames() );
if( !b )
{
    std::cerr << "Scanner failed" << std::endl;
    return 1;
}

gdcm::DiscriminateVolume dv;
dv.ProcessIntoVolume( s );
dv.Print( std::cout );

return 0;
}

```

## 29.36 DumbAnonymizer.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 This example shows how one can use the gdcm.Anonymizer in 'dumb' mode.
17 This class becomes really handy when one knows which particular tag to fill in.
18
19 Usage:
20
21 python DumbAnonymizer.py gdcmData/012345.002.050.dcm out.dcm
22
23 """
24
25 import gdcm
26
27 # http://www.oid-info.com/get/1.3.6.1.4.17434
28 THERALYS_ORG_ROOT = "1.3.6.1.4.17434"
29
30 tag_rules={
31     # Value
32     (0x0012,0x0010):("Value","MySponsorName"),
33     (0x0012,0x0020):("Value","MyProtocolID"),
34     (0x0012,0x0021):("Value","MyProtocolName"),
35     (0x0012,0x0062):("Value","YES"),
36     (0x0012,0x0063):("Value","MyDeidentificationMethod"),
37
38     # Method
39     (0x0002,0x0003):("Method","GenerateMSOPId"),
40     (0x0008,0x1155):("Method","GenerateMSOPId"),
41     (0x0008,0x0018):("Method","GenerateMSOPId"),
42     (0x0010,0x0010):("Method","GetSponsorInitials"),
43     (0x0010,0x0020):("Method","GetSponsorId"),

```

```

44 (0x0012,0x0030):("Method","GetSiteId"),
45 (0x0012,0x0031):("Method","GetSiteName"),
46 (0x0012,0x0040):("Method","GetSponsorId"),
47 (0x0012,0x0050):("Method","GetTPId"),
48 (0x0018,0x0022):("Method","KeepIfExist"),
49 (0x0018,0x1315):("Method","KeepIfExist"),
50 (0x0020,0x000d):("Method","GenerateStudyId"),
51 (0x0020,0x000e):("Method","GenerateSeriesId"),
52 (0x0020,0x1002):("Method","GetNumberOfFrames"),
53 (0x0020,0x0020):("Method","GetPatientOrientation"),
54 # Other:
55 (0x0012,0x0051):("Patient Field","Type Examen"),
56 (0x0018,0x1250):("Sequence Field","Receive Coil"),
57 (0x0018,0x0088):("Sequence Field","Spacing Between Slice"),
58 (0x0018,0x0095):("Sequence Field","Pixel Bandwidth"),
59 (0x0018,0x0082):("Sequence Field","Inversion Time"),
60 }
61
62 class MyAnon:
63     def __init__(self):
64         self.studyuid = None
65         self.seriesuid = None
66         generator = gdcm.UIDGenerator()
67         if not self.studyuid:
68             self.studyuid = generator.Generate()
69         if not self.seriesuid:
70             self.seriesuid = generator.Generate()
71     def GetSponsorInitials(self):
72         return "dummy^foobar"
73     def GenerateStudyId(self):
74         return self.studyuid
75     def GenerateSeriesId(self):
76         return self.seriesuid
77     #def GenerateMSOPIId(self):
78     def GenerateMSOPIId(self):
79         generator = gdcm.UIDGenerator()
80         return generator.Generate()
81     def GetSiteId(self):
82         return "MySiteId"
83     def GetSiteName(self):
84         return "MySiteName"
85     def GetSponsorId(self):
86         return "MySponsorId"
87     def GetTPId(self):
88         return "MyTP"
89
90 if __name__ == "__main__":
91     import sys
92     gdcm.FileMetaInformation.SetSourceApplicationEntityTitle
93     ( "DumbAnonymizer" )
94     gdcm.UIDGenerator.SetRoot( THERALYS_ORG_ROOT )
95
96 r = gdcm.Reader()
97 filename = sys.argv[1]
98 r.SetFileName( filename )
99 if not r.Read(): sys.exit(1)
100
101 obj = MyAnon()
102
103 w = gdcm.Writer()
104 ano = gdcm.Anonymizer()
105 ano.SetFile( r.GetFile() )
106 ano.RemoveGroupLength()
107 for tag,rule in tag_rules.items():
108     if rule[0] == 'Value':
109         print tag,rule
110         ano.Replace( gdcm.Tag( tag[0], tag[1] ), rule[1] )
111     elif rule[0] == 'Method':
112         print tag,rule
113         # result = locals()[rule[1]]()
114         methodname = rule[1]
115         if hasattr(obj, methodname):
116             _member = getattr(obj, methodname)
117             result = _member()
118             ano.Replace( gdcm.Tag( tag[0], tag[1] ), result )
119         else:
120             print "Problem with: ", methodname
121
122 outfilename = sys.argv[2]
123 w.SetFileName( outfilename )
124 w.SetFile( ano.GetFile() )

```

```
124     if not w.Write(): sys.exit(1)
```

## 29.37 DumpADAC.cxx

```
/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
/*
 * the goal of this example is to mimic the behavior of disp_img_header
 * see http://www.gmccorp-usa.com/IM/NM/GC/ADAC/SV/adactechtips/Released_01Q3.pdf
 */
#include "gdcmlReader.h"
#include "gdcmlPrivateTag.h"
#include "gdcmlAttribute.h"
#include "gdcmlImageWriter.h"

#include <iostream>
#include <fstream>
#include <vector>

#include <string.h>
#include <assert.h>
#include <stdint.h>

struct dict
{
    uint16_t key;
    const char *name;
};

dict Array[] = {
    { 0x01, "Patient name" },
    { 0x02, "Patient ID" },
    { 0x03, "Patient sex" },
    { 0x04, "Patient age" },
    { 0x05, "Patient height" },
    { 0x06, "Patient weight" },
    { 0x07, "Exam date" },
    { 0x08, "Dose admin. time" },
    { 0x09, "Unique exam key" },
    { 0x0a, "Exam procedure" },
    { 0x0b, "Referring physician" },
    { 0x0c, "Attending physician" },
    { 0x0d, "Imaging modality" },
    { 0x0e, "Hospital ID" },
    { 0x0f, "Histogram crv file" },
    { 0x10, "Acq. start time" },
    { 0x11, "Object data type" },
    { 0x12, "Image viewid" },
    { 0x13, "Imaging device name" },
    { 0x14, "Device serial number" },
    { 0x15, "Collimator" },
    { 0x16, "Software version" },
    { 0x17, "Radiopharmaceutical #1" },
    { 0x18, "Energy window #1 center" },
    { 0x19, "Radiopharmaceutical #2" },
    { 0x1a, "Energy window #1 width" },
    { 0x1b, "Isotope imaging mode" },
    { 0x1c, "Energy window #2 center" },
    { 0x1d, "Energy window #2 width" },
    { 0x1e, "Energy window #3 center" },
    { 0x1f, "Energy window #3 width" },
    { 0x20, "Energy window #4 center" },
    { 0x21, "Energy window #4 width" },
    { 0x22, "??Energy window #5 center" },
    { 0x23, "??Energy window #5 width" },
    { 0x24, "Patient orientation" },
```

```

{ 0x25, "Spatial resolution" },
{ 0x26, "Slice thickness" },
{ 0x27, "Image X dimension" },
{ 0x28, "Image Y dimension" },
{ 0x29, "Image Z dimension" },
{ 0x2a, "Image pixel width" },
{ 0x2b, "Uniformity corr. file" },
{ 0x2c, "Acquisition zoom factor" },
{ 0x2d, "Total counts in set" },
{ 0x2e, "Time / frame" },
{ 0x2f, "Total acq. time" },
{ 0x30, "Maximum pixel value" },
{ 0x31, "Minimum pixel value" },
{ 0x32, "R-R interval time" },
{ 0x33, "Percent of cycle imaged" },
{ 0x34, "# of cycles accepted" },
{ 0x35, "# of cycles rejected" },
{ 0x36, "Approximate ED frame" },
{ 0x37, "Approximate ES frame" },
{ 0x38, "Approximate EF" },
{ 0x39, "Starting angle" },
{ 0x3a, "Degrees of rotation" },
{ 0x3b, "Direction of rotation" },
{ 0x3c, "Cont. or step/shoot" },
{ 0x3d, "Lim recon start frame" },
{ 0x3e, "Upper window grey shade" },
{ 0x3f, "Lower lvl grey shade" },
{ 0x40, "Associated color map" },
{ 0x41, "Custom color map file" },
{ 0x42, "Manipulated image" },
{ 0x43, "Axis of rotation corr." },
{ 0x44, "Reorientation azimuth" },
{ 0x45, "Reorientation elevation" },
{ 0x46, "Filter type" },
{ 0x47, "Filter order" },
{ 0x48, "Filter cutoff frequency" },
{ 0x49, "Reconstruction type" },
{ 0x4a, "Attenuation coefficient" },
{ 0x4b, "Associated parent file" },
{ 0x4c, "Unique patient key" },
{ 0x52, "Normalization crv file" },
{ 0x53, "Unique object key" },
{ 0x54, "This phase of VFR is" },
{ 0x55, "True color value" },
{ 0x56, "# of sets of x,y,z grps" },
{ 0x57, "Scale factor of set" },
{ 0x6d, "Date of birth" },
{ 0x6e, "Directional orientation" },
{ 0x6f, "Number of VFR studies" },
{ 0x70, "R-R low tolerance" },
{ 0x71, "R-R high tolerance" },
{ 0x72, "Prog specific results:" },

{ 0x99, NULL }
};

void printname( int , int , uint16_t v )
{
    if( v == 0x1 )
    {
        std::cout << "DATABASE PARAMETERS" << std::endl;
        std::cout << "_____" << std::endl;
    }
    else if( v == 0x27 )
    {
        std::cout << "IMAGE PARAMETERS" << std::endl;
        std::cout << "_____" << std::endl;
    }
    else if( v == 0x13 )
    {
        std::cout << "EXTRA PARAMETERS" << std::endl;
        std::cout << "_____" << std::endl;
    }
    else if( v == 0x2e )
    {
        std::cout << "*** NOT CURRENTLY USED : " << std::endl;
    }
    static const unsigned int n = sizeof( Array ) / sizeof( *Array ) - 1;
    for( unsigned int i = 0; i < n; ++i )
    {
        if( v == Array[i].key )

```

```

        {
            std::cout << /*" " << std::dec << len << ", " << mult << " " << */ Array[i].name;
            std::cout << " : ";
            return;
        }
    }
    std::cout << /*"\t# " << std::dec << len << ", " << mult << */ std::hex << v << "\t: ";
}

uint16_t readint16(std::istream &is )
{
    uint16_t val;
    is.read( (char*)&val, sizeof( val ));
    return (uint16_t)((val>>8) | (val<<8));
}

uint32_t readint32(std::istream &is )
{
    uint32_t val;
    is.read( (char*)&val, sizeof( val ));
    val= ((val<<8)&0xFF00FF00) | ((val>>8)&0x00FF00FF);
    return (val>>16) | (val<<16);
}

float readfloat32(std::istream &is )
{
    union { uint32_t val; float f;} dual;
    dual.val = readint32(is);
    return dual.f;
}

struct el
{
    uint16_t v1;
    uint16_t v2;
    uint16_t v3;
    void read( std::istream & is )
    {
        v1 = readint16(is);
        v2 = readint16(is);
        v3 = readint16(is);
    }
    void print( std::ostream & os )
    {
        os << std::hex << v1 << "\t" << v2 << "\t" << v3 << std::endl;
    }
};

std::vector<el> Vel;

void readelement( std::istream & is )
{
    el e;
    e.read( is );
    Vel.push_back( e );
}

void printascii( uint16_t tag, const char *buffer, size_t len )
{
    std::ostream & os = std::cout;
    if( tag == 0x72 )
    {
        os << "\n ";
        for(size_t i = 0; i < len; ++i)
        {
            const char &c = buffer[i];
            if( c == 0x0 ) os << "!";
            else if( c == 0x0f ) os << " ";
            else if( c == 0x17 ) os << ":";
            else if( c == 0x14 ) os << ":";
            else if( c == 0x10 ) os << ":";
            else if( c == 0x16 ) os << ":";
            else if( c == 0x08 ) os << ":";
            else if( c == 0x0b ) os << ":";
            else if( c == 0x0e ) os << ":";
            else if( c == 0x07 ) os << ":";
            else os << c;
        }
        os << " ";
    }
    else

```

```

    {
        (void)len;
        os << "" << buffer << "";
    }
}

bool DumpADAC( std::istream & is )
{
    std::ostream &os = std::cout;

    char magic[6 + 1];
    magic[6] = 0;
    is.read( magic, 6);
    // std::cout << magic << " ";
    assert( strcmp( magic, "adac01" ) == 0 );
    int c = is.get();
    assert( c == 0 ); (void)c;
    c = is.get();
    assert( c == 'X' );

    uint16_t v;
    v = readint16(is);
    // std::cout << v << std::endl;
    assert( v == 512 ); (void)v; // ??

    int nel = 87;
    for (int i = 0; i <= nel; ++i )
    {
        readelement( is );
    }

    char buffer[512];
    for( int i = 0; i <= nel; ++i )
    {
        const el &e = Vel[i];
        int diff;
        if( i == nel )
        {
            diff = 2048 - e.v3;
            if( diff > 512 ) diff = 512;
        }
        else
        {
            const el &enext = Vel[i+1];
            diff = enext.v3 - e.v3;
        }
        is.seekg( e.v3, std::ios::beg );
        //std::cout << "(" << std::hex << std::setw( 2 ) << std::setfill( '0' ) << e.v1 << ")" " << std::hex <<
            std::setw( 3 ) << std::setfill( '0' ) << e.v2 << " ";
        printname( diff, 0, e.v1 );
        int mult = 1;
        if( e.v2 == 0 )
        {
            is.read( buffer, diff);
            buffer[ diff ] = 0;
            printascii( e.v1, buffer, diff);
        }
        else if( e.v2 == 0x100 )
        {
            mult = diff / 2;
            assert( diff == 2 * mult );
            for ( int ii = 0; ii < mult; ++ii )
            {
                if ( ii ) os << "\\ ";
                uint16_t val = readint16(is);
                os << "" << std::dec << val << "";
            }
        }
        else if( e.v2 == 0x200 )
        {
            assert( diff == 4 );
            uint32_t val = readint32(is);
            os << "" << std::dec << val << "";
        }
        else if( e.v2 == 0x300 )
        {
            assert( diff == 4 );
            float val = readfloat32(is);
            os << "" << std::dec << val << "";
        }
        else

```



```

        {
            assert( 0 );
        }
        os << std::endl;
    }
    return true;
}

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    // (0019,1061) UN (OB) 61\64\61\63\30          # 2048,1 Ver200 ADAC Pegasys Headers
    const gdcm::PrivateTag tver200adacpegasysheaders(0x0019,0x61,"ADAC_IMG");
    if( !ds.FindDataElement( tver200adacpegasysheaders ) ) return 1;
    const gdcm::DataElement& ver200adacpegasysheaders = ds.
        GetDataElement( tver200adacpegasysheaders );
    if ( ver200adacpegasysheaders.IsEmpty() ) return 1;
    const gdcm::ByteValue *bv = ver200adacpegasysheaders.
        GetByteValue();

    // (0019,1021) US 1                # 2,1 Ver200 Number of ADAC Headers
    // TODO

    // (0019,1041) IS [2048\221184 ] # 12,1-n Ver200 ADAC Header/Image Size
    if( bv->GetLength() != 2048 ) return 1;

    gdcm::Element<gdcm::VR::IS,gdcm::VM::VM2> el;
    const gdcm::PrivateTag tver200adacheaderimagesize(0x0019,0x41,"ADAC_IMG");
    if( !ds.FindDataElement( tver200adacheaderimagesize ) ) return 1;
    const gdcm::DataElement& ver200adacheaderimagesize = ds.
        GetDataElement( tver200adacheaderimagesize );
    el.SetFromDataElement( ver200adacheaderimagesize );
    if( el.GetValue(0) != 2048 ) return 1;

    std::stringstream is;
    std::string dup( bv->GetPointer(), bv->GetLength() );
    is.str( dup );
    bool b = DumpADAC( is );
    if( !b ) return 1;

    return 0;
}

```

## 29.38 DumpExamCard.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
/*

Try to extract contents of Philips RAW storage class:

(0002,0002) UI [1.2.840.10008.5.1.4.1.1.66]          # 26,1 Media Storage SOP Class UID
(0002,0003) UI [1.3.46.670589.11.17240.5.23.4.1.3012.2010032409482568018] # 56,1 Media Storage SOP
Instance UID

```

```

(0002,0010) UI [1.2.840.10008.1.2.1]          # 20,1 Transfer Syntax UID
(0002,0012) UI [1.3.46.670589.11.0.0.51.4.4.1] # 30,1 Implementation Class UID
(0002,0013) SH [MR DICOM 4.1]                 # 12,1 Implementation Version Name

* Everything done in this code is for the sole purpose of writing interoperable
* software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
* If you believe anything in this code violates any law or any of your rights,
* please contact us (gdcmm-developers@lists.sourceforge.net) so that we can
* find a solution.
*
* Everything you do with this code is at your own risk, since decompression
* algorithm was not written from specification documents.
*
* Special thanks to:
* Triplett, William T for bringing to your attention on this ExamCard stuff
*/
#include "gdcmmReader.h"
#include "gdcmmDataSet.h"
#include "gdcmmPrivateTag.h"
#include "gdcmmBase64.h"

#include <iomanip>

static bool compfn(const char *s1, const char *s2)
{
    return strcmp(s1,s2) < 0 ? true : false;
}

static const char *PDFStrings[] = { // Keep me ordered please
    "PDF_CONTROL_GEN_PARS",
    "PDF_CONTROL_PREP_PARS",
    "PDF_CONTROL_RECON_PARS",
    "PDF_CONTROL_SCAN_PARS",
    "PDF_EXAM_PARS",
    "PDF_HARDWARE_PARS",
    "PDF_PREP_PARS",
    "PDF_SPT_PARS",
};

static bool isvalidpdfstring( const char *pdfstring )
{
    assert( pdfstring );
    static const size_t n = sizeof( PDFStrings ) / sizeof( *PDFStrings );
    static const char **begin = PDFStrings;
    static const char **end = begin + n;
    return std::binary_search(begin, end, pdfstring, compfn);
}

typedef enum
{
    param_float = 0,
    param_integer,
    param_string,
    param_3, // ??
    param_enum,
} param_type;

static const char *gettypenamefromtype( int i)
{
    const char *ret = NULL;
    param_type e = (param_type)i;
    switch( e )
    {
        case param_float:
            ret = "float";
            break;
        case param_integer:
            ret = "int";
            break;
        case param_string:
            ret = "string";
            break;
        case param_3:
            ret = "??";
            break;
        case param_enum:
            ret = "enum";
            break;
    }
    assert( ret );
    return ret;
}

```

```

}

struct header
{
/*
 * TODO:
 * Looks as if we could read all int*, float* and string* at once...
 */
    int32_t v1; // offset to int pointer array ?
    uint16_t nints; // number of ints (max number?)
    uint16_t v3; // always 0 ?
    int32_t v4; // offset to float pointer array ?
    uint32_t nfloats;
    int32_t v6; // offset to string pointer array ?
    uint32_t nstrings;
    int32_t v8; // always 8 ??
    uint32_t numparams;
    uint32_t getnints() const { return nints; }
    uint32_t getnfloats() const { return nfloats; }
    uint32_t getnstrings() const { return nstrings; }
    uint32_t getnparams() const { return numparams; }
    void read( std::istream & is )
    {
        is.read( (char*)&v1, sizeof(v1));
        is.read( (char*)&nints, sizeof(nints));
        is.read( (char*)&v3, sizeof(v3));
        assert( v3 == 0 ); // looks like this is always 0
        is.read( (char*)&v4, sizeof(v4));
        is.read( (char*)&nfloats, sizeof(nfloats));
        is.read( (char*)&v6, sizeof(v6));
        is.read( (char*)&nstrings, sizeof(nstrings));
        is.read( (char*)&v8, sizeof(v8));
        assert( v8 == 8 );
        is.read( (char*)&numparams, sizeof(numparams));
    }
    void print( std::ostream & os )
    {
        os << v1 << ", ";
        os << nints << ", ";
        os << v3 << ", ";
        os << v4 << ", ";
        os << nfloats << ", ";
        os << v6 << ", ";
        os << nstrings << ", ";
        os << v8 << ", ";
        os << numparams << std::endl;
    }
};

struct param
{
    char name[32+1];
    int8_t boolean;
    int32_t type;
    uint32_t dim;
    uint32_t v4;
    /*int32_t*/ std::streamoff offset;
    param_type gettype() const { return (param_type)type; }
    uint32_t getdim() const { return dim; }
    void read( std::istream & is )
    {
        is.read( name, 32 + 1);
        //assert( name[32] == 0 ); // fails sometimes...
        // This is always the same issue the string can contains garbage from previous run,
        // we need to print only until the first \0 character:
        assert( strlen( name ) <= 32 ); // sigh
        is.read( (char*)&boolean, 1);
        assert( boolean == 0 || boolean == 1 ); // some kind of bool...
        is.read( (char*)&type, sizeof( type ) );
        assert( gettypenamefromtype( type ) );
        is.read( (char*)&dim, sizeof( dim ) );
        is.read( (char*)&v4, sizeof( v4 ) );
        //assert( v4 == 0 ); // always 0 ? sometimes not...
        const std::streamoff cur = is.tellg();
        is.read( (char*)&offset, sizeof( offset ) );
        offset += cur;
    }

    void print( std::ostream & os ) const
    {
        os << name << ", ";
    }
};

```

```

    os << (int)boolean << ", ";
    os << type << ", ";
    os << dim << ", ";
    os << v4 << ", ";
    os << offset << std::endl;
}

void printvalue( std::ostream & os, std::istream & is ) const
{
    is.seekg( offset );
    switch( type )
    {
        case param_float:
        {
            os.precision(2);
            os << std::fixed;
            for( uint32_t idx = 0; idx < dim; ++idx )
            {
                if( idx ) os << ", ";
                float v;
                is.read( (char*)&v, sizeof(v) );
                os << v; // what if the string contains \0 ?
            }
            break;
        }
        case param_integer:
        {
            for( uint32_t idx = 0; idx < dim; ++idx )
            {
                if( idx ) os << ", ";
                int32_t v;
                is.read( (char*)&v, sizeof(v) );
                os << v;
            }
            break;
        }
        case param_string:
        {
            std::string v;
            v.resize( dim );
            is.read( &v[0], dim );
            os << v;
        }
        break;
        case param_enum:
        {
            for( uint32_t idx = 0; idx < dim; ++idx )
            {
                if( idx ) os << ", ";
                int32_t v;
                is.read( (char*)&v, sizeof(v) );
                os << v;
            }
            break;
        }
    }
}

void printxml( std::ostream & os, std::istream & is ) const
{
    // <Attribute Name="CGEN_force_par_mode" Type="enum">0</Attribute>
    os << " <Attribute";
    os << " Name=\"" << name << "\"";
    os << " Type=\"" << gettypenamefromtype(type) << "\"";
    if( dim != 1 )
    {
        os << " ArraySize=\"" << dim << "\"";
    }
    os << ">";
    printvalue( os, is );
    os << "</Attribute>\n";
}

void printcsv( std::ostream & os, std::istream & is ) const
{
    os << std::setw(32) << std::left << name << ", ";
    os << std::setw(7) << std::right << gettypenamefromtype(type) << ", ";
    os << std::setw(4) << dim << ", ";
    os << " ";
    printvalue( os, is );
    os << ",\n";
}
};

```

```

static bool ProcessNested( gdcm::DataSet & ds )
{
    /*
    TODO:
    Looks like the real length of the blob is stored here:
(2005,1132) SQ                                     # u/1,1 ?
  (fffe,e000) na (Item with undefined length)
    (2005,0011) LO [Philips MR Imaging DD 002 ]      # 26,1 Private Creator
    (2005,1143) SL 3103                             # 4,1 ?

Wotsit ?
(2005,1132) SQ                                     # u/1,1 ?
  (fffe,e000) na (Item with undefined length)
    (2005,0011) LO [Philips MR Imaging DD 002 ]      # 26,1 Private Creator
    (2005,1147) CS [Y ]                             # 2,1 ?
*/
    bool ret = false;

    // (2005,1137) PN (LO) [PDF_CONTROL_GEN_PARS]      # 20,1 ?
    const gdcm::PrivateTag pt0(0x2005,0x37,"Philips MR Imaging DD 002");
    if( !ds.FindDataElement( pt0 ) ) return false;
    const gdcm::DataElement &de0 = ds.GetDataElement( pt0 );
    if( de0.IsEmpty() ) return false;
    const gdcm::ByteValue * bv0 = de0.GetByteValue();
    std::string s0( bv0->GetPointer() , bv0->GetLength() );

    // (2005,1139) LO [IEEE_PDF]                      # 8,1 ?
    const gdcm::PrivateTag pt1(0x2005,0x39,"Philips MR Imaging DD 002");
    if( !ds.FindDataElement( pt1 ) ) return false;
    const gdcm::DataElement &de1 = ds.GetDataElement( pt1 );

    const gdcm::PrivateTag pt(0x2005,0x44,"Philips MR Imaging DD 002");
    if( !ds.FindDataElement( pt ) ) return false;
    const gdcm::DataElement &de = ds.GetDataElement( pt );
    if( de.IsEmpty() ) return false;
    const gdcm::ByteValue * bv = de.GetByteValue();

    if( s0 == "ExamCardBlob" )
    {
        assert( de1.IsEmpty() );

        std::string fn = gdcm::LOComp::Trim( s0.c_str() ); // remove trailing space
        fn += ".xml";
        std::ofstream out( fn.c_str() );

        // remove trailing \0
        size_t len = strlen( bv->GetPointer() );
        out.write( bv->GetPointer() , len );
        out.close();

        // Extract binary64 thingy (this is a ugly hack, better use an XML parser)
        std::string dup( bv->GetPointer(), len );
        std::string::size_type pos1 = dup.find( "<ExamCardBlob>" );
        std::string::size_type pos2 = dup.find( "</ExamCardBlob>" );

        std::string b64( bv->GetPointer() + pos1 + 14, pos2 - (pos1 + 14) );

        // ugly hack to remove \r\n from input base64:
        std::string::iterator r_pos = std::remove(b64.begin(), b64.end(), '\r');
        b64.erase(r_pos, b64.end());
        std::string::iterator n_pos = std::remove(b64.begin(), b64.end(), '\n');
        b64.erase(n_pos, b64.end());
    }
    #if 0
        std::ofstream out2( "debug" );
        out2.write( b64.c_str(), b64.size() );
        out2.close();
    #endif

    const size_t dlen = gdcm::Base64::GetDecodeLength(b64.c_str(), b64.size() );

    std::string decoded;
    decoded.resize( dlen );
    gdcm::Base64::Decode( &decoded[0], decoded.size(), b64.c_str(), b64.size() );

    std::ofstream f64( "soap.xml" );
    f64.write( decoded.c_str(), decoded.size() );
    f64.close();

    ret = true;
}

```

```

    }
else
{
    if( del.IsEmpty() ) return false;
    const gdcm::ByteValue * bv1 = del.GetByteValue();
    std::string s1( bv1->GetPointer() , bv1->GetLength() );

    if( s1 == "IEEE_PDF" )
    {
        // std::cout << "Len= " << bv->GetLength() << std::endl;
    }
    #if 0
        std::string fn = gdcm::LOComp::Trim( s.c_str() ); // remove trailing space
        std::ofstream out( fn.c_str() );
        out.write( bv->GetPointer(), bv->GetLength() );
        out.close();
    #endif

    std::istreamstream is;
    std::string dup( bv->GetPointer(), bv->GetLength() );
    is.str( dup );

    header h;
    h.read( is );
    #if 0
        std::cout << s0.c_str() << std::endl;
        h.print( std::cout );
    #endif

    assert( is.tellg() == std::streampos(0x20) );
    is.seekg( 0x20 );

    std::vector< param > params;
    param p;
    for( uint32_t i = 0; i < h.getnparams(); ++i )
    {
        p.read( is );
        //p.print( std::cout );
        params.push_back( p );
    }

    std::string fn = gdcm::LOComp::Trim( s0.c_str() ); // remove trailing space
    assert( isvalidpdfstring( fn.c_str() ) );
    fn += ".csv";
    //fn += ".xml";
    std::ofstream csv( fn.c_str() );

    // let's do some bookkeeping:
    uint32_t nfloats = 0;
    uint32_t nints = 0;
    uint32_t nstrings = 0;
    for( std::vector<param>::const_iterator it = params.begin();
        it != params.end(); ++it )
    {
        param_type type = it->gettype();
        switch( type )
        {
            case param_float:
                nfloats += it->getdim();
                break;
            case param_integer:
                nints += it->getdim();
                break;
            case param_string:
                nstrings += it->getdim();
                break;
            default:
                ;
        }
    }
    #if 0
        std::cout << "Stats:" << std::endl;
        std::cout << "nfloats:" << nfloats << std::endl;
        std::cout << "nints:" << nints << std::endl;
        std::cout << "nstrings:" << nstrings << std::endl;
    #endif
    #endif

    assert( h.getnints() >= nints );
    assert( h.getnfloats() >= nfloats );
    assert( h.getnstrings() >= nstrings );

    for( uint32_t i = 0; i < h.getnparams(); ++i )
    {

```

```

        params[i].printcsv( csv, is );
        //params[i].printxml( csv, is );
    }
    csv.close();
    ret = true;
}
else if( s1 == "ASCII " )
{
    #if 0
        std::cerr << "ASCII is not handled" << std::endl;
        std::string fn = gdc::LOComp::Trim( s0.c_str() ); // remove trailing space
        fn += ".asc";
        std::ofstream out( fn.c_str() );
        out.write( bv->GetPointer() , bv->GetLength() );
        out.close();
    #endif

        std::string fn = gdc::LOComp::Trim( s0.c_str() ); // remove trailing space
        fn += ".sin";
        std::ofstream sin( fn.c_str() );

        const char *beg = bv->GetPointer();
        const char *end = beg + bv->GetLength();
        assert( *beg == 0 );
        const char *p = beg + 1; // skip first \0
        size_t prev = 0;
        for( ; p != end; ++p )
        {
            if( *p == 0 )
            {
                const char *s = beg + prev + 1;
                if( *s )
                {
                    sin << s << std::endl;
                }
                else
                {
                    sin << std::endl;
                }
                prev = p - beg;
            }
        }
        sin.close();

        ret = true;
    }
    else if( s1 == "BINARY" )
    {
        std::cerr << "BINARY is not handled" << std::endl;
        std::string fn = gdc::LOComp::Trim( s0.c_str() ); // remove trailing space
        fn += ".bin";
        std::ofstream out( fn.c_str() );
        //out.write( bv->GetPointer() + 512, bv->GetLength() - 512);
        out.write( bv->GetPointer() , bv->GetLength() );
        out.close();

        #if 0
            int array[ 128 ];
            memcpy( array, bv->GetPointer(), 512 );
            for( int i = 0; i < 14; ++i )
            {
                std::cout << array[i] << std::endl;
            }
        #endif

        ret = true;
    }
}
// else -> ret == false
assert( ret );

return ret;
}

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdc::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {

```

```

        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();
/*
(2005,1132) SQ                                     # u/1,1 ?
  (fffe,e000) na (Item with undefined length)
  (2005,0011) LO [Philips MR Imaging DD 002 ]      # 26,1 Private Creator
  (2005,1137) PN (LO) [PDF_CONTROL_GEN_PARS]      # 20,1 ?
  (2005,1138) PN (LO) (no value)                  # 0,1 ?
  (2005,1139) PN (LO) [IEEE_PDF]                  # 8,1 ?
  (2005,1140) PN (LO) (no value)                  # 0,1 ?
  (2005,1141) PN (LO) (no value)                  # 0,1 ?
  (2005,1143) SL 3103                             # 4,1 ?
  (2005,1144) OW
    66\05\00\00\3b\01\00\00\4a\0a\00\00\0e\00\00\00\7a\0a\00\00\95\01\00\00\08\00\00\00\1b\00\00\00\43\47\45\4e\5f\75\73\65\72\
    # 3104,1 ?
  (2005,1147) CS [Y ]                             # 2,1 ?
  (fffe,e00d)
*/
const gdcm::PrivateTag pt(0x2005,0x32,"Philips MR Imaging DD 002");
if( !ds.FindDataElement( pt ) ) return 1;
const gdcm::DataElement &de = ds.GetDataElement( pt );
if( de.IsEmpty() ) return 1;

gdcm::SequenceOfItems *sqi = de.GetValueAsSQ();
if ( !sqi ) return 1;
gdcm::SequenceOfItems::SizeType s = sqi->
  GetNumberOfItems();
for( gdcm::SequenceOfItems::SizeType i = 1; i <= s; ++i )
{
    gdcm::Item &item = sqi->GetItem(i);

    gdcm::DataSet &nestedds = item.GetNestedDataSet();

    if( !ProcessNested( nestedds ) ) return 1;
}

return 0;
}

```

## 29.39 DumpGEMSMovieGroup.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmImage.h"
#include "gdcmImageWriter.h"
#include "gdcmDataElement.h"
#include "gdcmPrivateTag.h"
#include "gdcmUIDGenerator.h"

#include <iostream>
#include <string>

#include <map>

bool PrintNameValuePairMapping( gdcm::SequenceOfItems *sqi_values,
gdcm::SequenceOfItems *sqi_names, std::string const & indent )
{
    using namespace gdcm;
    // prepare names mapping:
    typedef VRTToType<VR::UL>::Type UL;
    std::map< UL, std::string > names;
    assert( sqi_names );
    assert( sqi_values );

```



```

SequenceOfItems::SizeType s = sqi_names->
    GetNumberOfItems();
PrivateTag tindex(0x7fe1,0x71,"GEMS_Ultrasound_MovieGroup_001");
PrivateTag tname (0x7fe1,0x72,"GEMS_Ultrasound_MovieGroup_001");
// First sequence contains all possible names (this is a dict)
for( SequenceOfItems::SizeType i = 1; i <= s; ++i )
{
    const Item & item = sqi_names->GetItem( i );
    const DataSet & ds = item.GetNestedDataSet();
    if( !ds.FindDataElement( tindex )
        || !ds.FindDataElement( tname ) )
    {
        assert( 0 );
        return false;
    }
    const DataElement & index = ds.GetDataElement( tindex );
    const DataElement & name = ds.GetDataElement( tname );
    if( index.IsEmpty() || name.IsEmpty() )
    {
        assert( 0 );
        return false;
    }
    gdcmm::Element<VR::UL, VM::VM1> el1;
    el1.SetFromDataElement( index );

    gdcmm::Element<VR::LO, VM::VM1> el2;
    el2.SetFromDataElement( name );
    // std::cout << el1.GetValue() << " " << el2.GetValue() << std::endl;
    names.insert( std::make_pair( el1.GetValue(), el2.GetValue() ) );
}

SequenceOfItems::SizeType s2 = sqi_values->
    GetNumberOfItems();
assert( s2 <= s );
PrivateTag tindex2(0x7fe1,0x48,"GEMS_Ultrasound_MovieGroup_001");
for( SequenceOfItems::SizeType i = 1; i <= s2; ++i )
{
    const Item & item = sqi_values->GetItem( i );
    const DataSet & ds = item.GetNestedDataSet();
    if( !ds.FindDataElement( tindex2 ) )
    {
        assert( 0 );
        return false;
    }
    const DataElement & index2 = ds.GetDataElement( tindex2 );
    if( index2.IsEmpty() )
    {
        assert( 0 );
        return false;
    }
    gdcmm::Element<VR::FD, VM::VM1_2> el1;
    el1.SetFromDataElement( index2 );

    UL copy = (UL)el1.GetValue();
#ifdef 1
    std::cout << indent;
    std::cout << "( " << names[ copy ];
#endif
    // (7fe1,1052) FD 1560 # 8,1 ?
    // (7fe1,1057) LT [MscSkelSup] # 10,1 ?
    //PrivateTag tvalue(0x7fe1,0x52,"GEMS_Ultrasound_MovieGroup_001");
    PrivateTag tvalueint(0x7fe1,0x49,"GEMS_Ultrasound_MovieGroup_001"); // UL
    PrivateTag tvaluefloat1(0x7fe1,0x51,"GEMS_Ultrasound_MovieGroup_001"); // FL
    PrivateTag tvaluefloat(0x7fe1,0x52,"GEMS_Ultrasound_MovieGroup_001"); // FD
    PrivateTag tvalueul(0x7fe1,0x53,"GEMS_Ultrasound_MovieGroup_001"); // UL
    PrivateTag tvaluesl(0x7fe1,0x54,"GEMS_Ultrasound_MovieGroup_001"); // SL
    PrivateTag tvalueob(0x7fe1,0x55,"GEMS_Ultrasound_MovieGroup_001"); // OB
    PrivateTag tvaluetext(0x7fe1,0x57,"GEMS_Ultrasound_MovieGroup_001"); // LT
    PrivateTag tvaluefd(0x7fe1,0x77,"GEMS_Ultrasound_MovieGroup_001"); // FD / 1-N
    PrivateTag tvaluesl3(0x7fe1,0x79,"GEMS_Ultrasound_MovieGroup_001"); // SL / 1-N
    PrivateTag tvaluesl2(0x7fe1,0x86,"GEMS_Ultrasound_MovieGroup_001"); // SL ??
    PrivateTag tvaluefd1(0x7fe1,0x87,"GEMS_Ultrasound_MovieGroup_001"); // FD / 1-N
    PrivateTag tvaluefloat2(0x7fe1,0x88,"GEMS_Ultrasound_MovieGroup_001"); // FD ??
#ifdef 1
    std::cout << " ) = ";
#endif
    if( ds.FindDataElement( tvalueint ) )
    {
        const DataElement & value = ds.GetDataElement( tvalueint );
        gdcmm::Element<VR::UL,VM::VM1> el2;
        el2.SetFromDataElement( value );
    }
}

```

```

        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvaluefloat1 ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluefloat1 );
        gdcm::Element<VR::FL,VM::VM1> el2;
        el2.SetFromDataElement( value );
        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvaluefloat ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluefloat );
        gdcm::Element<VR::FD,VM::VM1> el2;
        el2.SetFromDataElement( value );
        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvalues1 ) )
    {
        const DataElement & value = ds.GetDataElement( tvalues1 );
        gdcm::Element<VR::SL,VM::VM1> el2;
        el2.SetFromDataElement( value );
        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvalueul ) )
    {
        const DataElement & value = ds.GetDataElement( tvalueul );
        gdcm::Element<VR::UL,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        assert( el2.GetLength() == 1 );
        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvalueob ) )
    {
        const DataElement & value = ds.GetDataElement( tvalueob );
        gdcm::Element<VR::SL,VM::VM1> el2;
        el2.SetFromDataElement( value );
        std::cout << el2.GetValue() << std::endl;
        std::cout << value << std::endl;
    }
    else if( ds.FindDataElement( tvaluetext ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluetext );
        gdcm::Element<VR::LT,VM::VM1> el2;
        el2.SetFromDataElement( value );
        std::cout << el2.GetValue() << std::endl;
    }
    else if( ds.FindDataElement( tvaluesl2 ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluesl2 );
        gdcm::Element<VR::SL,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        el2.Print( std::cout );
        assert( el2.GetLength() == 4 );
        std::cout << std::endl;
    }
    else if( ds.FindDataElement( tvaluesl3 ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluesl3 );
        gdcm::Element<VR::SL,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        el2.Print( std::cout );
        assert( el2.GetLength() == 4 );
        std::cout << std::endl;
    }
    else if( ds.FindDataElement( tvaluefd ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluefd );
        gdcm::Element<VR::FD,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        el2.Print( std::cout );
        assert( el2.GetLength() == 4 || el2.GetLength() == 3 || el2.GetLength() == 8 );
        std::cout << std::endl;
    }
    else if( ds.FindDataElement( tvaluefloat2 ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluefloat2 );
        gdcm::Element<VR::FD,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        el2.Print( std::cout );
        assert( el2.GetLength() == 2 );
        std::cout << std::endl;
    }

```

```

    }
    else if( ds.FindDataElement( tvaluefd1 ) )
    {
        const DataElement & value = ds.GetDataElement( tvaluefd1 );
        gdcmm::Element<VR::FD,VM::VM1_n> el2;
        el2.SetFromDataElement( value );
        el2.Print( std::cout );
        assert( el2.GetLength() == 4 );
        std::cout << std::endl;
    }
    else
    {
        std::cout << "(no value)" << std::endl;
        // std::cout << ds << std::endl;
        assert( ds.Size() == 2 );
    }
}
return true;
}

bool PrintNameValueMapping2( gdcmm::PrivateTag const & privtag, const
    gdcmm::DataSet & ds,
    gdcmm::SequenceOfItems *sqi_names, std::string const & indent )
{
    if( !ds.FindDataElement( privtag ) ) return 1;
    const gdcmm::DataElement& seq_values = ds.GetDataElement( privtag );
    gdcmm::SmartPointer<gdcmm::SequenceOfItems> sqi = seq_values.
        GetValueAssSQ();

    return PrintNameValueMapping( sqi, sqi_names, indent);
}

bool PrintNameValueMapping3( gdcmm::PrivateTag const & privtag1,
    gdcmm::PrivateTag const & privtag2, const gdcmm::DataSet & ds,
    gdcmm::SequenceOfItems *sqi_names, std::string const & indent )
{
    if( !ds.FindDataElement( privtag1 ) )
    {
        assert( 0 );
        return false;
    }
    const gdcmm::DataElement& values10name = ds.GetDataElement( privtag1 );
    gdcmm::Element<gdcmm::VR::LO,gdcmm::VM::VM1> el;
    el.SetFromDataElement( values10name );
    std::cout << std::endl;
    std::cout << " <" << el.GetValue().c_str() << ">" << std::endl;

    return PrintNameValueMapping2( privtag2, ds, sqi_names, indent);
}

bool print73( gdcmm::DataSet const & ds10, gdcmm::SequenceOfItems *sqi_dict,
    std::string const & indent )
{
    const gdcmm::PrivateTag tseq_values73(0x7fe1,0x73,"GEMS_Ultrasound_MovieGroup_001");
    if( !ds10.FindDataElement( tseq_values73 ) )
    {
        std::cout << indent << "No group 73" << std::endl;
        return false;
    }
    const gdcmm::DataElement& seq_values73 = ds10.GetDataElement( tseq_values73 );
    gdcmm::SmartPointer<gdcmm::SequenceOfItems> sqi_values73 =
        seq_values73.GetValueAssSQ();

    size_t ni3 = sqi_values73->GetNumberOfItems();
    for( size_t i3 = 1; i3 <= ni3; ++i3 )
    {
        gdcmm::Item &item_73 = sqi_values73->GetItem(i3);
        gdcmm::DataSet &ds73 = item_73.GetNestedDataSet();
        assert( ds73.Size() == 3 );

        const gdcmm::PrivateTag tseq_values74name(0x7fe1,0x74,"GEMS_Ultrasound_MovieGroup_001");
        const gdcmm::PrivateTag tseq_values75(0x7fe1,0x75,"GEMS_Ultrasound_MovieGroup_001");
        PrintNameValueMapping3( tseq_values74name, tseq_values75, ds73, sqi_dict, indent);
        std::cout << std::endl;
    }
    return true;
}

bool print36( gdcmm::DataSet const & ds10, gdcmm::SequenceOfItems *sqi_dict,
    std::string const & indent )

```

```

{
    const gdcm::PrivateTag tseq_values36(0x7fe1,0x36,"GEMS_Ultrasound_MovieGroup_001");
    if( !ds10.FindDataElement( tseq_values36 ) )
    {
        std::cout << indent << "No group 36" << std::endl;
        return false;
    }
    const gdcm::DataElement& seq_values36 = ds10.GetDataElement( tseq_values36
    );
    gdcm::SmartPointer<gdcm::SequenceOfItems> sqi_values36 =
        seq_values36.GetValueAsSQ();

    size_t ni3 = sqi_values36->GetNumberOfItems();
    assert( ni3 == 1 );
    for( size_t i3 = 1; i3 <= ni3; ++i3 )
    {
        gdcm::Item &item_36 = sqi_values36->GetItem(i3);
        gdcm::DataSet &ds36 = item_36.GetNestedDataSet();
        assert( ds36.Size() == 4 );

        // (7fe1,1037) UL 47 # 4,1 US MovieGroup Number of Frames
        // (7fe1,1043) OB 40\00\1c\c4\67\2f\0b\11\40 # 376,1 ?
        // (7fe1,1060) OB 4e\4e\49\4f\4e\47\46\43\2a # 4562714,1 US MovieGroup Image Data
        //
        const gdcm::PrivateTag timagedata(0x7fe1,0x60,"GEMS_Ultrasound_MovieGroup_001");
        assert( ds36.FindDataElement( timagedata ) );
        gdcm::DataElement const & imagedata = ds36.GetDataElement( timagedata );

        const gdcm::ByteValue * bv = imagedata.GetByteValue();
        assert( bv );
        static int c = 0;
        std::stringstream ss;
        ss << "/tmp/debug";
        ss << c++;
        std::ofstream os( ss.str().c_str(), std::ios::binary );
        os.write( bv->GetPointer(), bv->GetLength() );
        os.close();

        //const gdcm::PrivateTag tseq_values85(0x7fe1,0x85,"GEMS_Ultrasound_MovieGroup_001");
        //PrintNameValueMapping3( tseq_values84name, tseq_values85, ds83, sqi_dict, indent);
        //std::cout << std::endl;
    }
    return true;
}

bool print83( gdcm::DataSet const & ds10, gdcm::SequenceOfItems *sqi_dict
, std::string const & indent )
{
    const gdcm::PrivateTag tseq_values83(0x7fe1,0x83,"GEMS_Ultrasound_MovieGroup_001");
    if( !ds10.FindDataElement( tseq_values83 ) )
    {
        std::cout << indent << "No group 83" << std::endl;
        return false;
    }
    const gdcm::DataElement& seq_values83 = ds10.GetDataElement( tseq_values83
    );
    gdcm::SmartPointer<gdcm::SequenceOfItems> sqi_values83 =
        seq_values83.GetValueAsSQ();

    size_t ni3 = sqi_values83->GetNumberOfItems();
    for( size_t i3 = 1; i3 <= ni3; ++i3 )
    {
        gdcm::Item &item_83 = sqi_values83->GetItem(i3);
        gdcm::DataSet &ds83 = item_83.GetNestedDataSet();
        assert( ds83.Size() == 3 );

        const gdcm::PrivateTag tseq_values84name(0x7fe1,0x84,"GEMS_Ultrasound_MovieGroup_001");
        const gdcm::PrivateTag tseq_values85(0x7fe1,0x85,"GEMS_Ultrasound_MovieGroup_001");
        PrintNameValueMapping3( tseq_values84name, tseq_values85, ds83, sqi_dict, indent);
        std::cout << std::endl;
    }
    return true;
}

bool PrintNameValueMapping4( gdcm::PrivateTag const & privtag0, const
    gdcm::DataSet & subds, gdcm::PrivateTag const & privtag1,
    gdcm::PrivateTag const & privtag2,
    gdcm::SequenceOfItems *sqi_dict, std::string const & indent )
{
    (void)indent;
    if( !subds.FindDataElement( privtag0 ) )
    {

```

```

    assert( 0 );
    return 1;
}

const gdcm::DataElement& seq_values10 = subds.GetDataElement( privtag0 );
gdcm::SmartPointer<gdcm::SequenceOfItems> sqi_values10 =
    seq_values10.GetValueAsSQ();

size_t nil = sqi_values10->GetNumberOfItems();
// assert( nil == 1 );
for( size_t i1 = 1; i1 <= nil; ++i1 )
{
    gdcm::Item &item_10 = sqi_values10->GetItem(i1);
    gdcm::DataSet &ds10 = item_10.GetNestedDataSet();
    assert( ds10.Size() == 2 + 3 );
    // (7fe1,0010)
    // (7fe1,1012)
    // (7fe1,1018)
    // (7fe1,1020)
    // (7fe1,1083)

    PrintNameValueMapping3( privtag1, privtag2, ds10, sqi_dict, " " );
    std::cout << std::endl;

    const gdcm::PrivateTag tseq_values20(0x7fe1,0x20,"GEMS_Ultrasound_MovieGroup_001");
    if( !ds10.FindDataElement( tseq_values20 ) )
    {
        assert( 0 );
        return 1;
    }
    const gdcm::DataElement& seq_values20 = ds10.GetDataElement(
        tseq_values20 );
    gdcm::SmartPointer<gdcm::SequenceOfItems> sqi_values20 =
        seq_values20.GetValueAsSQ();

    size_t ni2 = sqi_values20->GetNumberOfItems();
    //assert( ni == 1 );
    for( size_t i2 = 1; i2 <= ni2; ++i2 )
    {
        gdcm::Item &item_20 = sqi_values20->GetItem(i2);
        gdcm::DataSet &ds20 = item_20.GetNestedDataSet();
        size_t count = ds20.Size(); (void)count;
        assert( ds20.Size() == 2 + 3 || ds20.Size() == 2 + 2 );
        // (7fe1,0010)
        // (7fe1,1024)
        // (7fe1,1026)
        // (7fe1,1036)
        // (7fe1,103a)
        // (7fe1,1083) (*)

        const gdcm::PrivateTag tseq_values20name(0x7fe1,0x24,"GEMS_Ultrasound_MovieGroup_001"
        );
        const gdcm::PrivateTag tseq_values26(0x7fe1,0x26,"GEMS_Ultrasound_MovieGroup_001");
        PrintNameValueMapping3( tseq_values20name, tseq_values26, ds20, sqi_dict, " " );
        std::cout << std::endl;

        print36(ds20, sqi_dict, " ");
        print83(ds20, sqi_dict, " ");
    }

    print83(ds10, sqi_dict, " ");
}
return true;
}

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    using namespace gdcm;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() ) return 1;

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();
    const PrivateTag tseq(0x7fe1,0x1,"GEMS_Ultrasound_MovieGroup_001");

    if( !ds.FindDataElement( tseq ) ) return 1;
    const DataElement& seq = ds.GetDataElement( tseq );

    SmartPointer<SequenceOfItems> sqi = seq.

```

```

    GetValueAsSQ();
    assert( sqi->GetNumberOfItems() == 1 );

    Item &item = sqi->GetItem(1);
    DataSet &subds = item.GetNestedDataSet();

    const PrivateTag tseq_dict(0x7fe1,0x70,"GEMS_Ultrasound_MovieGroup_001");
    if( !subds.FindDataElement( tseq_dict ) ) return 1;
    const DataElement& seq_dict = subds.GetDataElement( tseq_dict );
    SmartPointer<SequenceOfItems> sqi_dict = seq_dict.
        GetValueAsSQ();

    const PrivateTag tseq_values8(0x7fe1,0x8,"GEMS_Ultrasound_MovieGroup_001");
    if( !subds.FindDataElement( tseq_values8 ) ) return 1;
    const DataElement& seq_values8 = subds.GetDataElement( tseq_values8 );
    SmartPointer<SequenceOfItems> sqi_values8 = seq_values8.
        GetValueAsSQ();

    const PrivateTag tseq_values8name(0x7fe1,0x2,"GEMS_Ultrasound_MovieGroup_001");
    if( !subds.FindDataElement( tseq_values8name ) ) return 1;
    const DataElement& values8name = subds.GetDataElement( tseq_values8name );
{
    Element<VR::LO,VM::VM1> el;
    el.SetFromDataElement( values8name );
    std::cout << el.GetValue() << std::endl;
}
    size_t count = subds.Size(); (void)count;
    assert( subds.Size() == 3 + 2 + 1 || subds.Size() == 3 + 2 + 2);

// (7fe1,0010) # 30,1 Private Creator
// (7fe1,1002) # 8,1 US MovieGroup Value 0008 Name
// (7fe1,1003) # 4,1 ?
// (7fe1,1008) # 8140,1 US MovieGroup Value 0008 Sequence
// (7fe1,1010) # 1372196,1 ?
// (7fe1,1070) # 33684,1 US MovieGroup Dict
// (7fe1,1073) (*)
    PrintNameValueMapping( sqi_values8, sqi_dict, " ");

    const PrivateTag tseq_values10(0x7fe1,0x10,"GEMS_Ultrasound_MovieGroup_001");
    const PrivateTag tseq_values10name(0x7fe1,0x12,"GEMS_Ultrasound_MovieGroup_001");
    const PrivateTag tseq_values18(0x7fe1,0x18,"GEMS_Ultrasound_MovieGroup_001");
    PrintNameValueMapping4( tseq_values10, subds, tseq_values10name, tseq_values18, sqi_dict, " ");

    print73( subds, sqi_dict, " " );

#ifdef 0
    gdcm::DataSet::ConstIterator it = subds.Begin();
    for( ; it != subds.End(); ++it )
    {
        const gdcm::DataElement &de = *it;
        std::cout << de.GetTag() << std::endl;
    }
#endif

    return 0;
}

```

## 29.40 DumpImageHeaderInfo.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * Dump TOSHIBA MDW HEADER / Image Header Info
 */
#include "gdcmReader.h"
#include "gdcmPrivateTag.h"

```

Generated on Sat Sep 26 2015 07:03:31 for GDCM by Doxygen

```

    os << " " << val << std::endl;
    p += sizeof(val);
    char str2[17];
    memcpy( str2, p, 16 );
    str2[16] = 0;
    os << " " << str2 << std::endl;
}

#ifdef 0
std::ofstream out( str, std::ios::binary );
out.write( (char*)&magic, sizeof( magic ) );
out.write( (char*)&l, sizeof( l ) );
out.write( str, 16 );
out.write( &bytes[0], bytes.size() );
#endif
return is;
}

static bool DumpImageHeaderInfo( std::istream & is, size_t reflen )
{
    // TUSNONIMAGESTAM (5176)
    // TUSREMEASUREMEN (1352)
    // TUSBSINGLELAYOU (16)
    // TUSCLIPPARAMETE (104)

    element el;
    while( el.read( is ) )
    {
    }
    //size_t pos = is.tellg();
    //assert( pos == reflen );
    (void)reflen;

    return true;
}

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    const gdcm::PrivateTag timageheaderinfo(0x0029,0x10,"TOSHIBA MDW HEADER");
    if( !ds.FindDataElement( timageheaderinfo ) ) return 1;
    const gdcm::DataElement& imageheaderinfo = ds.GetDataElement(
        timageheaderinfo );
    if ( imageheaderinfo.IsEmpty() ) return 1;
    const gdcm::ByteValue * bv = imageheaderinfo.GetByteValue();

    std::stringstream is;
    std::string dup( bv->GetPointer(), bv->GetLength() );
    is.str( dup );
    bool b = DumpImageHeaderInfo( is, bv->GetLength() );
    if( !b ) return 1;

#ifdef 0
    const float d1 = 0.00416666668839752674; // 89 88 88 3B // 0x44c
    //const float d1 = 0.053231674455417881;
    const float d2 = 0.10828025639057159; // 0A C2 DD 3D // 0x1ac
    //const float d1 = 0.17869562069272813;
    //const unsigned int d2 = 4294967280;
    const float d3 = 0.10828025639057159; // 0A C2 DD 3D // 0x15c
    const int32_t d4 = 134;
    const uint32_t d5 = 1153476;
    std::ofstream t("/tmp/debug", std::ios::binary );
    //t.write( (char*)&d0, sizeof( d0 ) );
    t.write( (char*)&d1, sizeof( d1 ) );
    t.write( (char*)&d2, sizeof( d2 ) );
    t.write( (char*)&d3, sizeof( d3 ) );
    t.write( (char*)&d4, sizeof( d4 ) );
    t.write( (char*)&d5, sizeof( d5 ) );
    t.close();
#endif
}

```



```

    return 0;
}

```

## 29.41 DumpPhilipsECHO.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
#include "gdcmReader.h"
#include "gdcmDeflateStream.h"
#include "gdcm_zlib.h"

/*
 * This example extract the ZLIB compressed US image from a Philips private tag
 *
 * Everything done in this code is for the sole purpose of writing interoperable
 * software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
 * If you believe anything in this code violates any law or any of your rights,
 * please contact us (gdcm-developers@lists.sourceforge.net) so that we can
 * find a solution.
 *
 * Everything you do with this code is at your own risk, since decompression
 * algorithm was not written from specification documents.
 *
 * Usage:
 *
 * $ DumpPhilipsECHO private_us.dcm raw_us_img.raw
 * $ gdcming --sop-class-uid 1.2.840.10008.5.1.4.1.1.3.1 --size 608,427,88 raw_us_img.raw raw_us_img.dcm
 */

// header:
struct hframe
{
    uint32_t val0; // 800 increment ?
    uint16_t val1[2];
    uint16_t val2[2];
    uint32_t imgsize;

    bool operator==(const hframe &h) const
    {
        return val0 == h.val0 &&
            val1[0] == h.val1[0] &&
            val1[1] == h.val1[1] &&
            val2[0] == h.val2[0] &&
            val2[1] == h.val2[1] &&
            imgsize == h.imgsize;
    }
};

static bool ProcessDeflate( const char *outfilename, const int nslices, const
    int buf_size, const char *buf, const std::streampos len,
    const char *crdbuf, const size_t crclen )
{
    std::vector< hframe > crchheaders;
    crchheaders.reserve( nslices );
    {
        std::istringstream is;
        is.str( std::string( crdbuf, crclen ) );
        hframe header;
        for( int r = 0; r < nslices; ++r )
        {
            is.read( (char*)&header, sizeof( header ) );
        }
    }
    #if 0
        std::cout << header.val0
            << " " << header.val1[0]
            << " " << header.val1[1]

```

```

        << " " << header.val2[0]
        << " " << header.val2[1]
        << " " << header.imgsize << std::endl;
#endif
        crchheaders.push_back( header );
    }
}

std::istream is;
is.str( std::string( buf, len ) );

std::streamoff totalsize;
is.read( (char*)&totalsize, sizeof( totalsize ) );
assert( totalsize == len );

uint32_t nframes;
is.read( (char*)&nframes, sizeof( nframes ) );
assert( nframes == (uint32_t)nslices );

std::vector< std::streamoff > offsets;
offsets.reserve( nframes );
for( uint32_t frame = 0; frame < nframes ; ++frame )
{
    uint32_t offset;
    is.read( (char*)&offset, sizeof( offset ) );
    offsets.push_back( offset );
}

std::vector<char> outbuf;

const int size[2] = { 608, 427 }; // FIXME: where does it comes from ?
std::stringstream ss;
ss << outfilename;
ss << '_';
//ss << crchheaders[0].imgsize; // FIXME: Assume all header are identical !
ss << size[0];
ss << '_';
ss << size[1];
ss << '_';
ss << nframes;
ss << ".raw";
std::ofstream os( ss.str().c_str(), std::ios::binary );

assert( buf_size >= size[0] * size[1] );
outbuf.resize( buf_size );

hframe header;
//uint32_t prev = 0;
for( unsigned int r = 0; r < nframes; ++r )
{
    is.read( (char*)&header, sizeof( header ) );

    assert( header == crchheaders[r] );
    assert( header.val1[0] == 2000 );
    assert( header.val1[1] == 3 );
    assert( header.val2[0] == 1 );
    assert( header.val2[1] == 1280 );

    uLongf destLen = buf_size; // >= 608,427
    Bytef *dest = (Bytef*)&outbuf[0];
    assert( is.tellg() == offsets[r] + 16 );
    const Bytef *source = (Bytef*)buf + offsets[r] + 16;
    uLong sourceLen;
    if( r + 1 == nframes )
        sourceLen = totalsize - offsets[r] - 16;
    else
        sourceLen = offsets[r+1] - offsets[r] - 16;
    // FIXME: in-memory decompression:
    int ret = uncompress( dest, &destLen, source, sourceLen );
    assert( ret == Z_OK ); (void)ret;
    assert( destLen >= (uLongf)size[0] * size[1] ); // 16bytes padding ?
    assert( header.imgsize == (uint32_t)size[0] * size[1] );
    //os.write( &outbuf[0], outbuf.size() );
    os.write( &outbuf[0], size[0] * size[1] );

    // skip data:
    is.seekg( sourceLen, std::ios::cur );
}
os.close();
assert( is.tellg() == totalsize );

```

```

    return true;
}

static bool ProcessNone( const char *outfilename, const int nslices, const
    int buf_size, const char *buf, const std::streampos len,
    const char *crdbuf, const size_t crclen )
{
    std::vector< hframe > crcheaders;
    crcheaders.reserve( nslices );
    {
        std::istream is;
        is.str( std::string( crdbuf, crclen ) );
        hframe header;
        for( int r = 0; r < nslices; ++r )
        {
            is.read( (char*)&header, sizeof( header ) );
        }
    }
    #if 0
        std::cout << header.val0
            << " " << header.val1[0]
            << " " << header.val1[1]
            << " " << header.val2[0]
            << " " << header.val2[1]
            << " " << header.imgsize << std::endl;
    #endif
    crcheaders.push_back( header );
}

std::istream is;
is.str( std::string( buf, len ) );

std::streampos totalsize;
is.read( (char*)&totalsize, sizeof( totalsize ) );
assert( totalsize == len );

uint32_t nframes;
is.read( (char*)&nframes, sizeof( nframes ) );
assert( nframes == (uint32_t)nslices );

std::vector< uint32_t > offsets;
offsets.reserve( nframes );
for( uint32_t frame = 0; frame < nframes; ++frame )
{
    uint32_t offset;
    is.read( (char*)&offset, sizeof( offset ) );
    offsets.push_back( offset );
    //std::cout << offset << std::endl;
}

std::vector<char> outbuf;
// No idea how to present the data, I'll just append everything, and present it as 2D
std::stringstream ss;
ss << outfile;
ss << ' ';
ss << crcheaders[0].imgsize; // FIXME: Assume all header are identical !
ss << ' ';
ss << nframes;
ss << ".raw";
std::ofstream os( ss.str().c_str(), std::ios::binary );
outbuf.resize( buf_size ); // overallocated + 16
char *buffer = &outbuf[0];

hframe header;
for( unsigned int r = 0; r < nframes; ++r )
{
    is.read( (char*)&header, sizeof( header ) );
}
    #if 0
        std::cout << header.val0
            << " " << header.val1[0]
            << " " << header.val1[1]
            << " " << header.val2[0]
            << " " << header.val2[1]
            << " " << header.imgsize << std::endl;
    #endif
    assert( header == crcheaders[r] );

    is.read( buffer, buf_size - 16 );
    os.write( buffer, header.imgsize );
}
assert( is.tellg() == totalsize );
os.close();

```

```

    return true;
}

#ifdef NDEBUG
static const char * const UDM_USD_DATATYPE_STRINGS[] = {
    "UDM_USD_DATATYPE_DIN_2D_ECHO",
    "UDM_USD_DATATYPE_DIN_2D_ECHO_CONTRAST",
    "UDM_USD_DATATYPE_DIN_DOPPLER_CW",
    "UDM_USD_DATATYPE_DIN_DOPPLER_PW",
    "UDM_USD_DATATYPE_DIN_DOPPLER_PW_TDI",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_FLOW",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_PMI",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_CPA",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_TDI",
    "UDM_USD_DATATYPE_DIN_MMODE_ECHO",
    "UDM_USD_DATATYPE_DIN_MMODE_COLOR",
    "UDM_USD_DATATYPE_DIN_MMODE_COLOR_TDI",
    "UDM_USD_DATATYPE_DIN_PARAM_BLOCK",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_VELOCITY",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_POWER",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_VARIANCE",
    "UDM_USD_DATATYPE_DIN_DOPPLER_AUDIO",
    "UDM_USD_DATATYPE_DIN_DOPPLER_HIGHQ",
    "UDM_USD_DATATYPE_DIN_PHYSIO",
    "UDM_USD_DATATYPE_DIN_2D_COLOR_STRAIN",
    "UDM_USD_DATATYPE_DIN_COMPOSITE_RGB",
    "UDM_USD_DATATYPE_DIN_XFOV_REALTIME_GRAPHICS",
    "UDM_USD_DATATYPE_DIN_XFOV_MOSAIC",
    "UDM_USD_DATATYPE_DIN_COMPOSITE_R",
    "UDM_USD_DATATYPE_DIN_COMPOSITE_G",
    "UDM_USD_DATATYPE_DIN_COMPOSITE_B",
    "UDM_USD_DATATYPE_DIN_MMODE_COLOR_VELOCITY",
    "UDM_USD_DATATYPE_DIN_MMODE_COLOR_POWER",
    "UDM_USD_DATATYPE_DIN_MMODE_COLOR_VARIANCE",
    "UDM_USD_DATATYPE_DIN_2D_ELASTO",
};

static inline bool is_valid( const char * datatype_str )
{
    static const int n = sizeof( UDM_USD_DATATYPE_STRINGS ) / sizeof( *UDM_USD_DATATYPE_STRINGS );
    bool found = false;
    if( datatype_str )
    {
        for( int i = 0; !found && i < n; ++i )
        {
            found = strcmp( datatype_str, UDM_USD_DATATYPE_STRINGS[i] ) == 0;
        }
    }
    return found;
}
#endif

int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    using namespace gdcm;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() ) return 1;

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds1 = file.GetDataSet();

    const PrivateTag tseq1(0x200d,0x3cf8,"Philips US Imaging DD 045");
    if( !ds1.FindDataElement( tseq1 ) ) return 1;
    const DataElement& seq1 = ds1.GetDataElement( tseq1 );

    SmartPointer<SequenceOfItems> sq1 = seq1.
        GetValueAsSQ();
    assert( sq1->GetNumberOfItems() >= 1 );

    const size_t nitems = sq1->GetNumberOfItems();
    for( size_t item = 1; item < nitems; ++item )
    {
        Item &item1 = sq1->GetItem(item);
        DataSet &ds2 = item1.GetNestedDataSet();

        // (200d,300d) LO 28 UDM_USD_DATATYPE_DIN_2D_ECHO
        const PrivateTag tdatatype(0x200d,0x300d,"Philips US Imaging DD 033");

```

```

if( !ds2.FindDataElement( tdatatype ) ) return 1;
const DataElement& datatype = ds2.GetDataElement( tdatatype );
const ByteValue *bvdatatype = datatype.GetByteValue();
if( !bvdatatype ) return 1;

const PrivateTag tseq2(0x200d,0x3cf1,"Philips US Imaging DD 045");
if( !ds2.FindDataElement( tseq2 ) ) return 1;
const DataElement& seq2 = ds2.GetDataElement( tseq2 );

SmartPointer<SequenceOfItems> sqi2 = seq2.
    GetValueAsSQ();
assert( sqi2->GetNumberOfItems() >= 1 );

// FIXME: what if not in first Item ?
assert( sqi2->GetNumberOfItems() == 1 );
Item &item2 = sqi2->GetItem(1);
DataSet &ds3 = item2.GetNestedDataSet();

const PrivateTag tzlib(0x200d,0x3cfa,"Philips US Imaging DD 045");
if( !ds3.FindDataElement( tzlib ) ) return 1;
const DataElement& zlib = ds3.GetDataElement( tzlib );

const ByteValue *bv = zlib.GetByteValue();
if( !bv ) return 1;
if( bv->GetLength() != 4 ) return 1;

// (200d,3010) IS 2 88
const PrivateTag tnslices(0x200d,0x3010,"Philips US Imaging DD 033");
if( !ds3.FindDataElement( tnslices ) ) return 1;
const DataElement& nslices = ds3.GetDataElement( tnslices );
Element<VR::IS,VM::VM1> elnslices;
elnslices.SetFromDataElement( nslices );
const int nslicesref = elnslices.GetValue();
assert( nslicesref >= 0 );
// (200d,3011) IS 6 259648
const PrivateTag tzalloc(0x200d,0x3011,"Philips US Imaging DD 033");
if( !ds3.FindDataElement( tzalloc ) ) return 1;
const DataElement& zalloc = ds3.GetDataElement( tzalloc );
Element<VR::IS,VM::VM1> elzalloc;
elzalloc.SetFromDataElement( zalloc );
const int zallocref = elzalloc.GetValue();
assert( zallocref >= 0 );
// (200d,3021) IS 2 0
const PrivateTag tzzero(0x200d,0x3021,"Philips US Imaging DD 033");
if( !ds3.FindDataElement( tzzero ) ) return 1;
const DataElement& zero = ds3.GetDataElement( tzzero );
Element<VR::IS,VM::VM1> elzero;
elzero.SetFromDataElement( zero );
const int zerocref = elzero.GetValue();
assert( zerocref == 0 ); (void)zerocref;

// (200d,3cf3) OB
const PrivateTag tdeflate(0x200d,0x3cf3,"Philips US Imaging DD 045");
if( !ds3.FindDataElement( tdeflate ) ) return 1;
const DataElement& deflate = ds3.GetDataElement( tdeflate );
const ByteValue *bv2 = deflate.GetByteValue();

// (200d,3cfb) OB
const PrivateTag tcrc(0x200d,0x3cfb,"Philips US Imaging DD 045");
if( !ds3.FindDataElement( tcrc ) ) return 1;
const DataElement& crc = ds3.GetDataElement( tcrc );
const ByteValue *bv3 = crc.GetByteValue();

std::string outfile = std::string( bvdatatype->GetPointer(), bvdatatype->
    GetLength() );
outfile = LOComp::Trim( outfile.c_str() );
const char *outfilename = outfile.c_str();
assert( is_valid(outfilename) );
if( bv2 )
{
    assert( bv3 );
    assert( zallocref > 0 );
    assert( nslicesref > 0 );
    std::cout << ds2 << std::endl;

    if( strncmp(bv->GetPointer(), "ZLib", 4) == 0 )
    {
        if( !ProcessDeflate( outfile, nslicesref, zallocref, bv2->GetPointer(),
            std::streampos(bv2->GetLength()), bv3->GetPointer(), bv3->
            GetLength() ) )
        {

```

```

        return 1;
    }
}
else if( strcmp(bv->GetPointer(), "None", 4) == 0 )
{
    if( !ProcessNone( outfilename, nslicesref, zallocref, bv2->GetPointer(),
        std::streampos(bv2->GetLength()), bv3->GetPointer(), bv3->
        GetLength() ) )
    {
        return 1;
    }
}
else
{
    std::string str( bv->GetPointer(), bv->GetLength() );
    std::cerr << "Unhandled: " << str << std::endl;
    return 1;
}
}
}

return 0;
}

```

## 29.42 DumpToSQLITE3.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * Ref:
 * http://massmail.spl.harvard.edu/public-archives/slicer-devel/2010/004408.html
 *
 * Implementation details:
 * http://www.sqlite.org/c3ref/bind_blob.html
 * http://www.adp-gmbh.ch/sqlite/bind_insert.html
 */
#include "gdcmScanner.h"
#include "gdcmDirectory.h"
#include "gdcmTag.h"
#include "gdcmTrace.h"

#include "sqlite3.h"

#include <stdio.h>
#include <time.h>

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        return 1;
    }
    time_t time_start = time(0);

    gdcm::Trace::SetDebug( false );
    gdcm::Trace::SetWarning( false );
    const char *inputdirectory = argv[1];

    gdcm::Directory d;
    unsigned int nfiles = d.Load( inputdirectory, true);

    gdcm::Scanner s;
    using gdcm::Tag;
    s.AddTag( Tag(0x20,0xd) ); // Study Instance UID
    s.AddTag( Tag(0x20,0xe) ); // Series Instance UID

```

```

bool b0 = s.Scan( d.GetFileNames() );
if( !b0 ) return 1;
time_t time_scanner = time(0);

std::cout << "Finished loading data from : " << nfiles << " files" << std::endl;

// MappingType const &mappings = s.GetMappings();

sqlite3* db;
sqlite3_open("./dicom.db", &db);

if(db == 0)
{
    std::cerr << "Could not open database." << std::endl;
    return 1;
}

const char sql_stmt[] = "create table browser (seriesuid, studyuid)";
int ret;

char *errmsg;
ret = sqlite3_exec(db, sql_stmt, 0, 0, &errmsg);

if(ret != SQLITE_OK)
{
    printf("Error in statement: %s [%s].\n", sql_stmt, errmsg);
    return 1;
}
using gdcmm::Directory;
using gdcmm::Scanner;
const Directory::FileNamesType& files = d.GetFileNames();
Directory::FileNamesType::const_iterator file = files.begin();

sqlite3_stmt *stmt;
if ( sqlite3_prepare(
    db,
    "insert into browser values (?,?)", // stmt
    -1, // If than zero, then stmt is read up to the first nul terminator
    &stmt,
    0 // Pointer to unused portion of stmt
)
!= SQLITE_OK)
{
    printf("\nCould not prepare statement.");
    return 1;
}
//printf("\nThe statement has %d wildcards\n", sqlite3_bind_parameter_count(stmt));
for(; file != files.end(); ++file)
{
    const char *filename = file->c_str();
    bool b = s.IsKey(filename);
    if( b )
    {
        const Scanner::TagToValue &mapping = s.GetMapping(filename);
        Scanner::TagToValue::const_iterator it = mapping.begin();

        sqlite3_reset(stmt);

        for( int index = 1; it != mapping.end(); ++it, ++index)
        {
            //const Tag &tag = it->first;
            const char *value = it->second;

            if (sqlite3_bind_text (
                stmt,
                index, // Index of wildcard
                value,
                (int)strlen(value), // length of text
                SQLITE_STATIC // SQLite assumes that the information is in static
            )
            != SQLITE_OK)
            {
                printf("\nCould not bind int.\n");
                return 1;
            }
        }
        if (sqlite3_step(stmt) != SQLITE_DONE)
        {
            printf("\nCould not step (execute) stmt.\n");
            return 1;
        }
    }
}

```

```

    }
}

sqlite3_close(db);

time_t time_sqlite = time(0);

std::cout << "Time to scan DICOM files: " << (time_scanner - time_start) << std::endl;
std::cout << "Time to build SQLITE3: " << (time_sqlite - time_scanner) << std::endl;

return 0;
}

```

## 29.43 DuplicatePCDE.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmlReader.h"
#include "gdcmlWriter.h"
#include "gdcmlItem.h"
#include "gdcmlImageReader.h"
#include "gdcmlSequenceOfItems.h"
#include "gdcmlFile.h"
#include "gdcmlTag.h"
/*
Usage:
DuplicatePCDE gdcmlData/D_CLUNIE_CT1_J2KI.dcm out.dcm

aka:
medical.nema.org/medical/dicom/DataSets/WG04/IMAGES/J2KI/CT1_J2KI

See:
gdcmlConformanceTests/CT1_J2KI_DuplicatePCDE.dcm

Original thread can be found at:

http://groups.google.com/group/comp.protocols.dicom/browse_thread/thread/82f28c4db28963af

Question:
1.
There is no restriction for a specific Private Creator Data Element
(PCDE) to be unique within the same group, right ?
Decoders of Private Data would have to handle the case where a PCDE
would be repeated and should NOT stop on the first instance of a
particular PCDE, right ?

Eg. when searching for the tag associated with
(0x0029,0x0010,"SIEMENS CSA HEADER") in the following (pseudo)
dataset:

(0029,0010) LO [SIEMENS CSA HEADER] # 18, 1
PrivateCreator
(0029,0011) LO [SIEMENS MEDCOM HEADER] # 22, 1
PrivateCreator
(0029,0012) LO [SIEMENS MEDCOM HEADER2] # 22, 1
PrivateCreator
(0029,0013) LO [SIEMENS CSA HEADER] # 18, 1
PrivateCreator
(0029,1008) CS [IMAGE NUM 4] # 12, 1
CSAImageHeaderType
(0029,1009) LO [20050723] # 8, 1
CSAImageHeaderVersion
(0029,1010) OB 53\56\31\30\04\03\02\01\38\00\00\00\4d
\00\00\00\45\63\68\6f\4c\69... # 6788, 1 CSAImageHeaderInfo

```



```

(0029,1018) CS [MR] # 2, 1
CSASeriesHeaderType
(0029,1019) LO [20050723] # 8, 1
CSASeriesHeaderVersion
(0029,1020) OB 53\56\31\30\04\03\02\01\2c\00\00\00\4d
\00\00\00\55\73\65\64\50\61... # 51520, 1 CSASeriesHeaderInfo
(0029,1131) LO [4.0.163088300] # 14, 1
PMTFInformation1
(0029,1132) UL 32768 # 4, 1
PMTFInformation2
(0029,1133) UL 0 # 4, 1
PMTFInformation3
(0029,1134) CS [DB TO DICOM] # 12, 1
PMTFInformation4
(0029,1260) ?? 63\6f\6d\20 # 4, 1
Unknown Tag & Data
(0029,1310) OB 53\56\31\30\04\03\02\01\38\00\00\00\4d
\00\00\00\45\63\68\6f\4c\69... # 6788, 1 CSAImageHeaderInfo

```

one should return two instances, correct ?

Answer:

I would say that this is covered in principle by the PS 3.5 7.1  
 "The Data Elements ... shall occur at most once in a Data Set"  
 rule, since the data element is defined by the tuple  
 (private creator,gggg,ee) where xxee is the element  
 number and xx is arbitrary and has no inherent meaning and  
 does not serve to disambiguate the data element.

E.g.:

```

(0019,0030) Private Creator ID = "Smith"
...
(0019,0032) Private Creator ID = "Smith"
...
(0019,3015) Fractal Index = "32"
...
(0019,3215) Fractal Index = "32"

```

would be illegal because even though they are assigned different  
 (completely arbitrary) blocks, with the same group, element  
 number and private creator, (0019,3015) and (0019,3215) are the  
 "same" data element.

\*/

```

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();

    // Let's get all private element from group 0x9:
    /*
(0009,0010) LO [GEMS_IDEN_01] # 12,1 Private Creator
(0009,1001) LO [GE_GENESIS_FF ] # 14,1 Full fidelity
(0009,1002) SH [CT01] # 4,1 Suite id
(0009,1004) SH [HiSpeed CT/i] # 12,1 Product id
(0009,1027) SL 862399669 # 4,1 Image actual date
(0009,1030) SH (no value) # 0,1 Service id
(0009,1031) SH (no value) # 0,1 Mobile location number
(0009,10e6) SH [05] # 2,1 Genesis Version - now
(0009,10e7) UL 973283917 # 4,1 Exam Record checksum
(0009,10e9) SL 862399669 # 4,1 Actual series data time stamp
*/
    gdcm::Tag start(0x0009,0x0);
    // Create a temporary duplicate dataset, since we cannot insert data element as we go over them (std::set
    // would reorganize itself as we go over it ...)

```

```

gdcmm::DataSet dup;
gdcmm::Tag new_private(0x0009,0x0);
while (start.GetGroup() == 0x9 )
{
    const gdcmm::DataElement& de = ds.FindNextDataElement(start);
    const gdcmm::Tag &t = de.GetTag();
    if( t.IsPrivateCreator() )
    {
        std::cout << t << std::endl;
        // Ok let's duplicate into the next available attribute:
        gdcmm::DataElement duplicate = de;
        duplicate.GetTag().SetElement( (uint16_t)(t.GetElement() + 1) );
        dup.Insert( duplicate );
        new_private = duplicate.GetTag();
    }
    else if( t.IsPrivate() && !t.IsPrivateCreator() )
    {
        //std::cout << de << std::endl;
        std::string owner = ds.GetPrivateCreator( de.GetTag() );
        //std::cout << owner << std::endl;
        gdcmm::DataElement duplicate = de;
        duplicate.GetTag().SetPrivateCreator( new_private );
        if( const gdcmm::ByteValue *bv = duplicate.GetByteValue() )
        {
            // Warning: when doing : duplicate = de, only the pointer to the ByteValue is passed
            // (to avoid large memory duplicate). We need to explicitly duplicate the bytevalue ourselves:
            gdcmm::ByteValue *dupbv = new gdcmm::ByteValue( bv->GetPointer(),
                bv->GetLength() );
            // Let's recognize the duplicated ASCII-type elements:
            if( duplicate.GetVR() & gdcmm::VR::VRASCII )
                dupbv->Fill( 'X' );
            duplicate.SetValue( *dupbv );
        }
        dup.Insert( duplicate );
    }
    start = t;
    // move to next possible 'public' element
    start.SetElement( (uint16_t)(start.GetElement() + 1) );
}

gdcmm::DataSet::ConstIterator it = dup.Begin();
for( ; it != dup.End(); ++it )
{
    ds.Insert( *it );
}

gdcmm::Writer w;
w.SetFile( file );
w.SetFileName( outfilename );
if ( !w.Write() )
{
    return 1;
}

return 0;
}

```

## 29.44 ELSCINT1WaveToText.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmmReader.h"
#include "gdcmmPrivateTag.h"

/*
 * This example shows how to read a Wave Information tag from ELSCINT1

```

```

* The wave information is stored in Tag (01e1,18,ELSCINT1) hidden in a
* Secondary Capture Image Storage (usually a 'N' Symbol is shown)
*
* Everything done in this code is for the sole purpose of writing interoperable
* software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
* If you believe anything in this code violates any law or any of your rights,
* please contact us (gdcmm-developers@lists.sourceforge.net) so that we can
* find a solution.
*
* Everything you do with this code is at your own risk, since decompression
* algorithm was not written from specification documents.
*
* Special thanks to:
* Gauthier Bouilhol
*/

template <typename T>
bool dumpargs(std::ostream & os, T c1, T c2, T c3, T c4, T c5, T c6, T c7, T c8)
{
    static const char sep = '\t';
    os << c1 << sep << c2 << sep << c3 << sep << c4 << sep << c5 << sep << c6 << sep << c7 << sep << c8;
    os << std::endl;
    return true;
}

bool wave2stream( std::ostream &text_file, const char *in, size_t len )
{
    short * buffer = (short*)in;
    size_t length = len / sizeof( short );
    text_file << "COMPLETE_WAVE" << '\t' << "MASK" << '\t' << "AQUSITION_PROFIL" << '\t' << "
    END-INHALE" << '\t' << "END-EXHALE" << '\t' << "AQUSITION_WAVE" << '\t' << "WAVE_STATISTICS" << '\t' << "MASK"
    << std::endl;
    for (size_t i=0;i<length-76;i+=2)
    {
        if ( i < 74 )
        {
            if (buffer[i+75] == 0)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            if (buffer[i+75] == 16384)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            if (buffer[i+75] == 256)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << buffer[i+74] << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            if (buffer[i+75] == -32768)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << buffer[i+74] << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            if (buffer[i+75] == -16384)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << buffer[i+74] << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            if (buffer[i+75] == -32512)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << buffer[i+74] << '\t' << buffer[i+74] << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
        }
        else
        {
            if (buffer[i+75] == 0)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << " " << '\t' << " " << std::endl;
            if (buffer[i+75] == 16384)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << " " << '\t' << " " << std::endl;
            if (buffer[i+75] == 256)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 0 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << " " << std::endl;
            if (buffer[i+75] == -32768)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << " " << std::endl;
            if (buffer[i+75] == -16384)
                text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " " << '\t' << buffer[i] << '\t' << buffer[i+1] << std::endl;
            buffer[i+74] << '\t' << " " << '\t' << " " << '\t' << " " << '\t' << " " << std::endl;
        }
    }
}

```

```

        buffer[i+74] << '\t' << " " << '\t' << buffer[i+74] << '\t' << " " << '\t' << " "
        << std::endl;
        if (buffer[i+75] == -32512)
            text_file << buffer[i+74] << '\t' << buffer[i+75] << '\t' << 1 << '\t' << " "
            << '\t' << buffer[i+74] << '\t' << buffer[i+74] << '\t' << " " << '\t' << " "
            << std::endl;
    }
}

return true;
}

int main(int argc, char *argv [])
{
    if( argc < 3 ) return 1;
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    const gdcm::PrivateTag twave(0x01e1,0x18,"ELSCINT1");
    if( !ds.FindDataElement( twave ) ) return 1;
    const gdcm::DataElement& wave = ds.GetDataElement( twave );
    if ( wave.IsEmpty() ) return 1;
    const gdcm::ByteValue * bv = wave.GetByteValue();
    assert( bv );

    std::ofstream os( outfile, std::ios::binary );
    // Dump that to a CSV file:
    wave2stream( os, bv->GetPointer(), bv->GetLength() );
    os.close();

    return 0;
}

```

## 29.45 EncapsulateFileInRawData.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmAnonymizer.h"
#include "gdcmWriter.h"
#include "gdcmUIDGenerator.h"
#include "gdcmFile.h"
#include "gdcmTag.h"
#include "gdcmSystem.h"

#include "magic.h" // libmagic, API to file command line tool

/*
 * Let say you want to encapsulate a file type that is not defined in DICOM (exe, zip, png)
 * PNG is a bad example, unless it contains transparency (which has been deprecated).
 * It will take care of dispatching each chunk to an appropriate data item (pretty much like
 * WaveformData)
 *
 * Usage:
 * ./EncapsulateFileInRawData large_input_file.exe large_input_file.dcm
 */

// TODO:
// $ file -bi /tmp/gdcm-2.1.0.pdf

```

```

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " inputfile output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    if( !gdcmm::System::FileExists( filename ) ) return 1;

    size_t s = gdcmm::System::FileSize(filename);
    if( !s ) return 1;

    magic_t cookie = magic_open(MAGIC_NONE);
    const char * file_type = magic_file(cookie, filename);
    if( !file_type ) return 1;
    magic_close(cookie);

    gdcmm::Writer w;
    gdcmm::File &file = w.GetFile();
    //gdcmm::DataSet &ds = file.GetDataSet();
    //w.SetCheckFileMetaInformation( true );
    w.SetFileName( outfile );

    file.GetHeader().SetDataSetTransferSyntax(
        gdcmm::TransferSyntax::ImplicitVRLittleEndian );

    gdcmm::Anonymizer anon;
    anon.SetFile( file );

    gdcmm::MediaStorage ms = gdcmm::MediaStorage::RawDataStorage
        ;

    gdcmm::UIDGenerator gen;
    anon.Replace( gdcmm::Tag(0x0008,0x16), ms.GetString() );
    std::cout << ms.GetString() << std::endl;
    anon.Replace( gdcmm::Tag(0x0008,0x18), gen.Generate() );

    if( !w.Write() )
    {
        std::cerr << "Could not write: " << outfile << std::endl;
        return 1;
    }

    return 0;
}

```

## 29.46 ExtractEncapsulatedFile.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * This example shows how one from C# context can extract a binary blob
 * and write out as a file.
 * This example is meant for pdf encapsulated file, but can be adapted for other type
 * of binary blob.
 *
 * DICOM file is:
 * ...
 * (0042,0010) ST (no value available) # 0, 0 DocumentTitle
 * (0042,0011) OB 25\50\44\46\2d\31\2e\32\20\0d\25\2e\3cf\d3\20\0d\31\30\20\30\20... # 40718, 1
 * EncapsulatedDocument

```

```

* (0042,0012) LO [application/pdf] # 16, 1 MIMETimeTypeOfEncapsulatedDocument
* ...
*
* Usage:
* $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
* $ mono bin/ExtractEncapsulatedFile.exe some_pdf_encapsulated.dcm
*/
using System;
using gdcm;

public class ExtractEncapsulatedFile
{
    public static int Main(string[] args)
    {
        string file = args[0];
        Reader reader = new Reader();
        reader.SetFileName( file );
        bool ret = reader.Read();
        if( !ret )
        {
            return 1;
        }

        File f = reader.GetFile();
        DataSet ds = f.GetDataSet();
        Tag tencapsulated_stream = new Tag(0x0042,0x0011); // Encapsulated Document
        if( !ds.FindDataElement( tencapsulated_stream ) )
        {
            return 1;
        }
        // else
        DataElement de = ds.GetDataElement( tencapsulated_stream );
        ByteValue bv = de.GetByteValue();
        uint len = bv.GetLength();
        byte[] encapsulated_stream = new byte[len];
        bv.GetBuffer( encapsulated_stream, len );

        // Write out the decompressed bytes
        //System.Console.WriteLine(image.toString());
        using (System.IO.Stream stream =
            System.IO.File.Open(@"tmp/dd.pdf",
                System.IO.FileMode.Create))
        {
            System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
            writer.Write( encapsulated_stream );
        }

        return 0;
    }
}

```

## 29.47 ExtractEncryptedContent.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
#include "gdcmReader.h"

#include <fstream>

/*

openssl smime -encrypt -binary -aes256 -in outputfile.dcm -inform DER -out outputfile.der -outform DER ../
trunk/Testing/Source/Data/certificate.pem

openssl smime -decrypt -binary -in out.der -inform DER -out outputfile.dcm -outform DER -inkey ../trunk/

```

```

        Testing/Source/Data/privatekey.pem ../trunk/Testing/Source/Data/certificate.pem

    */

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.der" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();

    const gdcm::DataElement &EncryptedAttributesSequence = ds.
        GetDataElement( gdcm::Tag( 0x0400,0x0500 ) );

    gdcm::SequenceOfItems *sqi = EncryptedAttributesSequence.
        GetValueAsSQ();

    if ( !sqi || sqi->GetNumberOfItems() != 1 ) return 1;

    gdcm::Item &item = sqi->GetItem(1);

    gdcm::DataSet &nesteddds = item.GetNestedDataSet();

    if( ! nesteddds.FindDataElement( gdcm::Tag( 0x0400,0x0520 ) ) ) return 1;

    const gdcm::DataElement &EncryptedContent = nesteddds.
        GetDataElement( gdcm::Tag( 0x0400,0x0520 ) );

    const gdcm::ByteValue *bv = EncryptedContent.GetByteValue();

    std::ofstream of( outfile, std::ios::binary );
    of.write( bv->GetPointer(), bv->GetLength() );
    of.close();

    return 0;
}

```

## 29.48 ExtractIconFromFile.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example shows how to either retrieve an Icon if present somewhere
 * in the file, or else generate one.
 */
#include "gdcmImageReader.h"
#include "gdcmPNMCodec.h"
#include "gdcmIconImageFilter.h"
#include "gdcmIconImageGenerator.h"

bool WriteIconAsPNM(const char* filename, const gdcm::IconImage& icon)
{
    gdcm::PNMCodec pnm;

```

```

    pnm.SetDimensions( icon.GetDimensions() );
    pnm.SetPixelFormat( icon.GetPixelFormat() );
    pnm.SetPhotometricInterpretation( icon.
        GetPhotometricInterpretation() );
    pnm.SetLUT( icon.GetLUT() );
    const gdcm::DataElement& in = icon.GetDataElement();
    bool b = pnm.Write( filename, in );
    assert( b );
    return b;
}

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcm::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read (or not image): " << filename << std::endl;
        return 1;
    }

    gdcm::IconImageFilter iif;
    iif.SetFile( reader.GetFile() );
    bool b = iif.Extract();

    if( b )
    {
        const gdcm::IconImage &icon = iif.GetIconImage(0);
        icon.Print( std::cout );

        if( !icon.GetTransferSyntax().IsEncapsulated() )
        {
            // Let's write out this icon as PNM file
            WriteIconAsPNM("icon.ppm", icon);
        }
        else if( icon.GetTransferSyntax() ==
            gdcm::TransferSyntax::JPEGBaselineProcess1
            || icon.GetTransferSyntax() ==
            gdcm::TransferSyntax::JPEGExtendedProcess2_4
        )
        {
            const gdcm::DataElement& in = icon.GetDataElement();
            const gdcm::ByteValue *bv = in.GetByteValue();
            assert( bv );
            std::ofstream out( "icon.jpg", std::ios::binary );
            out.write( bv->GetPointer(), bv->GetLength() );
            out.close();
        }
    }
    else
    {
        assert( iif.GetNumberOfIconImages() == 0 );
        std::cerr << "No Icon Found anywhere in file" << std::endl;

        const gdcm::Image &img = reader.GetImage();
        gdcm::IconImageGenerator iig;
        iig.AutoPixelMinMax(true);
        iig.SetPixmap( img );
        const unsigned int idims[2] = { 64, 64 };
        iig.SetOutputDimensions( idims );
        //iig.SetPixelMinMax(60, 868);
        if( !iig.Generate() ) return 1;
        const gdcm::IconImage &icon = iig.GetIconImage();
        WriteIconAsPNM("icon.ppm", icon);
    }

    return 0;
}

```

## 29.49 ExtractImageRegion.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

```



Copyright (c) 2006-2011 Mathieu Malaterre  
 All rights reserved.  
 See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even  
 the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR  
 PURPOSE. See the above copyright notice for more information.

```

=====*/

/*
 * This small code shows how to use the gdcm.ImageRegionReader API
 * In this example we are taking each frame by frame and dump them to
 * /tmp/frame.raw.
 *
 * Usage:
 * $ bin/ExtractImageRegion.exe input.dcm
 *
 * Example:
 * $ bin/ExtractImageRegion.exe gdcmData/012345.002.050.dcm
 * $ md5sum /tmp/frame.raw
 * d594a5e2fde12f32b6633ca859b4d4a6 /tmp/frame.raw
 * $ gdcminfo --md5sum gdcmData/012345.002.050.dcm
 * [...]
 * md5sum: d594a5e2fde12f32b6633ca859b4d4a6
 */
using System;
using gdcm;

public class ExtractImageRegion
{
    public static int Main(string[] args)
    {
        string filename = args[0];

        uint file_size = gdcm.PosixEmulation.FileSize(filename);

        // instantiate the reader:
        gdcm.ImageRegionReader reader = new gdcm.
            ImageRegionReader();
        reader.SetFileName( filename );

        // pull DICOM info:
        if (!reader.ReadInformation()) return 1;

        // store current offset:
        uint cur_pos = reader.GetStreamCurrentPosition();

        uint remaining = file_size - cur_pos;

        Console.WriteLine("Remaining bytes to read (Pixel Data): " + remaining.ToString() );

        // Get file infos
        gdcm.File f = reader.GetFile();

        // get some info about image
        UIntArrayType dims = ImageHelper.GetDimensionsValue(f);
        PixelFormat pf = ImageHelper.GetPixelFormatValue( f);
        int pixelsize = pf.GetPixelSize();
        PhotometricInterpretation pi = ImageHelper.GetPhotometricInterpretationValue(f);
        Console.WriteLine( pi.ToString() );

        // buffer to get the pixels
        byte[] buffer = new byte[ dims[0] * dims[1] * pixelsize ];

        // define a simple box region.
        BoxRegion box = new BoxRegion();
        for (uint z = 0; z < dims[2]; z++)
        {
            // Define that I want the image 0, full size (dimx x dimy pixels)
            // and do that for each z:
            box.SetDomain(0, dims[0] - 1, 0, dims[1] - 1, z, z);
            //System.Console.WriteLine( box.ToString() );
            reader.SetRegion( box );

            // reader will try to load the uncompressed image region into buffer.
            // the call returns an error when buffer.Length is too small. For instance
            // one can call:
            // uint buf_len = reader.ComputeBufferLength(); // take into account pixel size
            // to get the exact size of minimum buffer
            if (reader.ReadIntoBuffer(buffer, (uint)buffer.Length))

```

```

        {
            using (System.IO.Stream stream =
                System.IO.File.Open(@"tmp/frame.raw",
                    System.IO.FileMode.Create))
            {
                System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
                writer.Write(buffer);
            }
        }
    else
    {
        throw new Exception("can't read pixels error");
    }
}

return 0;
}
}

```

## 29.50 ExtractImageRegion.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

/*
 * This small code shows how to use the gdcm.ImageRegionReader API
 * In this example we are taking each frame by frame and dump them to
 * /tmp/frame.raw.
 *
 * Usage:
 * $ LD_LIBRARY_PATH=. CLASSPATH=gdcm.jar:. java ExtractImageRegion input.dcm
 */
import gdcm.*;
import java.io.FileOutputStream;

public class ExtractImageRegion
{
    public static void main(String[] args) throws Exception
    {
        String filename = args[0];

        // instantiate the reader:
        ImageRegionReader reader = new ImageRegionReader();
        reader.SetFileName( filename );

        // pull DICOM info:
        if (!reader.ReadInformation()) return;
        // Get file infos
        File f = reader.GetFile();

        // get some info about image
        UIntArrayType dims = ImageHelper.GetDimensionsValue(f);
        PixelFormat pf = ImageHelper.GetPixelFormatValue(f);
        int pixelsize = pf.GetPixelSize();

        // buffer to get the pixels
        long buffer_length = dims.get(0) * dims.get(1) * pixelsize;
        byte[] buffer = new byte[ (int)buffer_length ];

        // define a simple box region.
        BoxRegion box = new BoxRegion();
        for (int z = 0; z < dims.get(2); z++)
        {
            // Define that I want the image 0, full size (dimx x dimy pixels)
            // and do that for each z:
            box.SetDomain(0, dims.get(0) - 1, 0, dims.get(1) - 1, z, z);

```

```

//System.Console.WriteLine( box.toString() );
reader.SetRegion( box );

// reader will try to load the uncompressed image region into buffer.
// the call returns an error when buffer.Length is too small. For instance
// one can call:
// long buf_len = reader.ComputeBufferLength(); // take into account pixel size
// to get the exact size of minimum buffer
if (reader.ReadIntoBuffer(buffer, buffer_length))
{
    FileOutputSteam fos = new FileOutputSteam("/tmp/frame.raw");
    fos.write(buffer);
    fos.close();
}
else
{
    throw new Exception("can't read pixels error");
}
}
}

```

## 29.51 ExtractImageRegionWithLUT.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * This small code shows how to use the gdcml.ImageRegionReader API
 * In this example we are taking each frame by frame and dump them to
 * /tmp/frame.raw.
 * Furthermore we are applying the LUT on this image.
 * Special care should be taken in case the image is not PALETTE COLOR
 *
 * Usage:
 * $ bin/ExtractImageRegionWithLUT.exe input.dcm
 *
 * Example:
 * $ bin/ExtractImageRegionWithLUT.exe gdcmlData/rle16l00.dcm
 * $ md5sum /tmp/frame_rgb.raw
 * 73bf61325fdb6e2830244a2b7b0c4ae2 /tmp/frame_rgb.raw
 * $ gdcmlimg --depth 16 --spp 3 --size 600,430 /tmp/frame_rgb.raw rgb.dcm
 * $ gdcmlviewer rgb.dcm
 */
using System;
using gdcml;

public class ExtractImageRegion
{
    public static int Main(string[] args)
    {
        string filename = args[0];

        // instantiate the reader:
        gdcml.ImageRegionReader reader = new gdcml.
            ImageRegionReader();
        reader.SetFileName( filename );

        // pull DICOM info:
        if (!reader.ReadInformation()) return 1;
        // Get file infos
        gdcml.File f = reader.GetFile();

        gdcml.LookupTable lut = reader.GetImage().GetLUT();
    }
}

```

```

// get some info about image
UIntArrayType dims = ImageHelper.GetDimensionsValue(f);
PixelFormat pf = ImageHelper.GetPixelFormatValue (f);
int pixelsize = pf.GetPixelSize();

// buffer to get the pixels
byte[] buffer = new byte[ dims[0] * dims[1] * pixelsize ];

// output buffer for the RGB decoded image:
byte[] buffer2 = new byte[ dims[0] * dims[1] * pixelsize * 3 ];

// define a simple box region.
BoxRegion box = new BoxRegion();
for (uint z = 0; z < dims[2]; z++)
{
    // Define that I want the image 0, full size (dimx x dimy pixels)
    // and do that for each z:
    box.SetDomain(0, dims[0] - 1, 0, dims[1] - 1, z, z);
    //System.Console.WriteLine( box.ToString() );
    reader.SetRegion( box );

    // reader will try to load the uncompressed image region into buffer.
    // the call returns an error when buffer.Length is too small. For instance
    // one can call:
    // uint buf_len = reader.ComputeBufferLength(); // take into account pixel size
    // to get the exact size of minimum buffer
    if (reader.ReadIntoBuffer(buffer, (uint)buffer.Length))
    {
        if( !lut.Decode( buffer2, (uint)buffer2.Length, buffer, (uint)buffer.Length ) )
        {
            throw new Exception("can't decode");
        }

        using (System.IO.Stream stream =
            System.IO.File.Open(@"tmp/frame_rgb.raw",
                System.IO.FileMode.Create))
        {
            System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
            writer.Write(buffer2);
        }
    }
    else
    {
        throw new Exception("can't read pixels error");
    }
}

return 0;
}
}

```

## 29.52 Extracting\_All\_Resolution.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
// This work was realised during the GSOC 2011 by Manoj Alwani

#include <fstream>
#include <openjpeg.h>
#include <stdint.h>
#include <string.h>
#include <assert.h>
#include <gdcm_j2k.h>
#include <gdcm_jp2.h>
#include <iostream>
#include <cstring>

```

```

#include <stdio.h>
#include <string.h>
#include <stdlib.h>
#include <math.h>
#include "gdcmImageReader.h"
#include "gdcmSequenceOfItems.h"
#include "gdcmSystem.h"
#include <fstream>

#include "gdcmMediaStorage.h"
#include "gdcmWriter.h"
#include "gdcmItem.h"
#include "gdcmImageReader.h"
#include "gdcmAttribute.h"
#include "gdcmFile.h"
#include "gdcmTag.h"
#include "gdcmTransferSyntax.h"
#include "gdcmUIDGenerator.h"
#include "gdcmAnonymizer.h"
#include "gdcmStreamImageWriter.h"
#include "gdcmImageHelper.h"
#include "gdcmTrace.h"

void error_callback(const char *msg, void *) {
    (void)msg;
}

void warning_callback(const char *msg, void *) {
    (void)msg;
}

void info_callback(const char *msg, void *) {
    (void)msg;
}

bool Write_Resolution(gdcm::StreamImageWriter & theStreamWriter, const char *
    filename, int res, std::ostream& of, int flag, gdcm::SequenceOfItems *sq, int
    No_Of_Resolutions)
{
    std::ifstream is;
    is.open( filename, std::ios::binary );
    opj_dparameters_t parameters; /* decompression parameters */
    opj_event_mgr_t event_mgr; /* event manager */
    opj_dinfo_t* dinfo; /* handle to a decompressor */
    opj_cio_t *cio;
    opj_image_t *image = NULL;
    // FIXME: Do some stupid work:
    is.seekg( 0, std::ios::end);
    std::streampos buf_size = is.tellg();
    char *dummy_buffer = new char[(unsigned int)buf_size];
    is.seekg(0, std::ios::beg);
    is.read( dummy_buffer, buf_size);
    unsigned char *src = (unsigned char*)dummy_buffer;
    uint32_t file_length = (uint32_t)buf_size; // 32bits truncation should be ok since DICOM cannot have
        larger than 2Gb image

    /* configure the event callbacks (not required) */
    memset(&event_mgr, 0, sizeof(opj_event_mgr_t));
    event_mgr.error_handler = error_callback;
    event_mgr.warning_handler = warning_callback;
    event_mgr.info_handler = info_callback;

    /* set decoding parameters to default values */
    opj_set_default_decoder_parameters(&parameters);

    // default blindly copied
    parameters.cp_layer=0;
    parameters.cp_reduce= res;
    // parameters.decode_format=-1;
    // parameters.cod_format=-1;

    const char jp2magic[] = "\x00\x00\x00\x0C\x6A\x50\x20\x20\x0D\x0A\x87\x0A";
    if( memcmp( src, jp2magic, sizeof(jp2magic) ) == 0 )
    {
        /* JPEG-2000 compressed image data ... sigh */
        // gdcmData/ELSCINT1_JP2vsJ2K.dcm
        // gdcmData/MAROTECH_CT_JP2Lossy.dcm
        //gdcmWarningMacro( "J2K start like JPEG-2000 compressed image data instead of codestream" );
        parameters.decode_format = 1; //JP2_CFMT;
    }
}

```

```

    //assert(parameters.decode_format == JP2_CFMT);
}
else
{
    /* JPEG-2000 codestream */
    //parameters.decode_format = J2K_CFMT;
    //assert(parameters.decode_format == J2K_CFMT);
    assert( 0 );
}
parameters.cod_format = 11; // PGX_DFMT;
//assert(parameters.cod_format == PGX_DFMT);

/* get a decoder handle */
dinfo = opj_create_decompress(CODEC_JP2);

/* catch events using our callbacks and give a local context */
opj_set_event_mgr((opj_common_ptr)dinfo, &event_mgr, NULL);

/* setup the decoder decoding parameters using user parameters */
opj_setup_decoder(dinfo, &parameters);

/* open a byte stream */
cio = opj_cio_open((opj_common_ptr)dinfo, src, file_length);

/* decode the stream and fill the image structure */
image = opj_decode(dinfo, cio);
if(!image) {
    opj_destroy_decompress(dinfo);
    opj_cio_close(cio);
    //gdcErrorMacro( "opj_decode failed" );
    return 1;
}

    opj_cp_t * cp = ((opj_jp2_t*)dinfo->jp2_handle)->j2k->cp;
    opj_tcp_t *tcp = &cp->tcps[0];
    opj_tccp_t *tccp = &tcp->tccps[0];
    /* std::cout << "\n No of Cols In Image" << image->x1;
    std::cout << "\n No of Rows In Image" << image->y1;
    std::cout << "\n No of Components in Image" << image->numcomps;
    std::cout << "\n No of Resolutions"<< tccp->numresolutions << "\n";
*/
    opj_j2k_t* j2k = NULL;
    opj_jp2_t* jp2 = NULL;
    jp2 = (opj_jp2_t*)dinfo->jp2_handle;
    int reversible = jp2->j2k->cp->tcps->tccps->qmfbid;
    //std::cout << reversible;
    int compno = 0;
    opj_image_comp_t *comp = &image->comps[compno];
    int Dimensions[2];
    Dimensions[0]= comp->w;
    Dimensions[1] = comp->h;
    opj_cio_close(cio);
    unsigned long len = Dimensions[0]*Dimensions[1] * image->numcomps;
    //std::cout << "\nTest" <<image->comps[0].factor;
    char *raw = new char[len];
    for (unsigned int compno = 0; compno < (unsigned int)image->numcomps; compno++)
    {
        opj_image_comp_t *comp = &image->comps[compno];

        int w = image->comps[compno].w;
        int h = image->comps[compno].h;
        uint8_t *data8 = (uint8_t*)raw + compno;
        for (int i = 0; i < w * h ; i++)
        {
            int v = image->comps[compno].data[i];
            *data8 = (uint8_t)v;
            data8 += image->numcomps;
        }
    }

    gdc::Writer w;
    gdc::File &file = w.GetFile();
    gdc::DataSet &ds = file.GetDataSet();

    file.GetHeader().SetDataSetTransferSyntax(
        gdc::TransferSyntax::ExplicitVRLittleEndian );

    gdc::UIDGenerator uid;
    gdc::DataElement de( gdc::Tag(0x8,0x18) ); // SOP Instance UID
    de.SetVR( gdc::VR::UI );
    const char *u = uid.Generate();

```

```

de.SetByteValue( u, strlen(u) );
ds.Insert( de );

gdcm::DataElement de1( gdcm::Tag(0x8,0x16) );
de1.SetVR( gdcm::VR::UI );
gdcm::MediaStorage ms( gdcm::MediaStorage::CTImageStorage
);
de1.SetByteValue( ms.GetString(), strlen(ms.GetString()) );
ds.Insert( de1 );

const char mystr[] = "MONOCHROME2 ";
gdcm::DataElement de2( gdcm::Tag(0x28,0x04) );
//de.SetTag( gdcm::Tag(0x28,0x04) );
de2.SetVR( gdcm::VR::CS );
de2.SetByteValue(mystr, strlen(mystr));
ds.Insert( de2 );

gdcm::Attribute<0x0028,0x0010> row = {image->comps[0].w};
//row.SetValue(512);
ds.Insert( row.GetAsDataElement() );
// w.SetCheckFileMetaInformation( true );
gdcm::Attribute<0x0028,0x0011> col = {image->comps[0].h};
ds.Insert( col.GetAsDataElement() );
gdcm::Attribute<0x0028,0x0008> Number_Of_Frames = {1};
ds.Insert( Number_Of_Frames.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0100> at = {8};
ds.Insert( at.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0002> at1 = {image->numcomps};
ds.Insert( at1.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0101> at2 = {8};
ds.Insert( at2.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0102> at3 = {7};
ds.Insert( at3.GetAsDataElement() );

if (flag == 1)
{
    for (int i=0; i < No_Of_Resolutions; i++)
    {
        int a = 1;
        int b = 1;

        while(a!=((No_Of_Resolutions)-i))
        {
            b = b*2;
            a = a+1;
        }
        uint16_t row = (image->y1)/b;
        uint16_t col = (image->x1)/b;
        //std::cout << row;
        gdcm::Element<gdcm::VR::IS, gdcm::VM::VM1> el2;
        el2.SetValue(i+1);
        gdcm::DataElement rfn = el2.GetAsDataElement(); //ulr --> upper
            left row
        rfn.SetTag( gdcm::Tag(0x0008,0x1160) );

        gdcm::Element<gdcm::VR::US, gdcm::VM::VM2> el;
        el.SetValue(1,0);
        el.SetValue(1,1);
        gdcm::DataElement ulr = el.GetAsDataElement(); //ulr --> upper
            left col/row
        ulr.SetTag( gdcm::Tag(0x0048,0x0201) );

        gdcm::Element<gdcm::VR::US, gdcm::VM::VM2> el1;
        el1.SetValue(col,0);
        el1.SetValue(row,1);
        gdcm::DataElement brr = el1.GetAsDataElement();
        brr.SetTag( gdcm::Tag(0x0048,0x0202) ); //brr --> bottom right col/row
        gdcm::Item it;
        gdcm::DataSet &nds = it.GetNestedDataSet();
        nds.Insert( rfn );
        nds.Insert(ulr);
        nds.Insert(brr);
    }
}

```

```

    sq->AddItem(it);
}

gdcm::Writer w1;
gdcm::File &file1 = w1.GetFile();
gdcm::DataSet &ds1 = file1.GetDataSet();
file1.GetHeader().SetDataSetTransferSyntax(
    gdcm::TransferSyntax::ExplicitVRLittleEndian );

gdcm::UIDGenerator uid1;
gdcm::DataElement dea( gdcm::Tag(0x8,0x18) ); // SOP Instance UID
dea.SetVR( gdcm::VR::UI );
const char *u1 = uid1.Generate();
dea.SetByteValue( u1, strlen(u1) );
ds1.Insert( dea );

gdcm::DataElement deb( gdcm::Tag(0x8,0x16) );
deb.SetVR( gdcm::VR::UI );
gdcm::MediaStorage ms1(
    gdcm::MediaStorage::VLWholeSlideMicroscopyImageStorage
);
deb.SetByteValue( ms1.GetString(), strlen(ms1.GetString()) );
ds1.Insert( deb );

const char mystr1[] = "MONOCHROME2 ";
gdcm::DataElement dec( gdcm::Tag(0x28,0x04) );
//de.SetTag(gdcm::Tag(0x28,0x04));
dec.SetVR( gdcm::VR::CS );
dec.SetByteValue(mystr, strlen(mystr));
ds1.Insert( dec );

gdcm::Attribute<0x0028,0x0010> row1 = {image->y1};
//row.SetValue(512);
ds1.Insert( row1.GetAsDataElement() );
// w.SetCheckFileMetaInformation( true );
gdcm::Attribute<0x0028,0x0011> col1 = {image->x1};
ds1.Insert( col1.GetAsDataElement() );
gdcm::Attribute<0x0028,0x0008> Number_Of_Frames1 = {tccp->numresolutions};
ds1.Insert( Number_Of_Frames1.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0100> ata = {8};
ds1.Insert( ata.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0002> atb = {image->numcomps};
ds1.Insert( atb.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0101> atc = {8};
ds1.Insert( atc.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0102> atd = {7};
ds1.Insert( atd.GetAsDataElement() );

theStreamWriter.SetFile(file1);

gdcm::DataElement des( gdcm::Tag(0x0048,0x0200) );
des.SetVR(gdcm::VR::SQ);
//des.SetVR(gdcm::VM::VM1);
des.SetValue(*sq);
des.SetVLToUndefined();

ds1.Insert( des );

if (!theStreamWriter.WriteImageInformation()){
    std::cerr << "unable to write image information" << std::endl;
    return 1; //the CanWrite function should prevent getting here, else,
    //that's a test failure
}

}

theStreamWriter.SetFile(file);

if (!theStreamWriter.CanWriteFile()){
    delete [] raw;
    std::cout << "Not able to write";
    return 0; //this means that the file was unwritable, period.
    //very similar to a ReadImageInformation failure
}
else
    std::cout<<"\nable to read";

```



```

// Important to write here
std::vector<unsigned int> extent = gdcm::ImageHelper::GetDimensionsValue
    (file);

unsigned short xmax = extent[0];
unsigned short ymax = extent[1];
unsigned short theChunkSize = 4;
unsigned short ychunk = extent[1]/theChunkSize; //go in chunk sizes of theChunkSize
unsigned short zmax = extent[2];
std::cout << "\n" << xmax << "\n" << ymax << "\n" << zmax << "\n" << image->numcomps << "\n";

if (xmax == 0 || ymax == 0)
{
    std::cerr << "Image has no size, unable to write zero-sized image." << std::endl;
    return 0;
}

int z, y, nexty;
unsigned long prevLen = 0; //when going through the char buffer, make sure to grab
//the bytes sequentially. So, store how far you got in the buffer with each iteration.
for (z = 0; z < zmax; ++z){
    for (y = 0; y < ymax; y += ychunk){
        nexty = y + ychunk;
        if (nexty > ymax) nexty = ymax;
        theStreamWriter.DefinePixelExtent(0, xmax, y, nexty, z, z+1);
        unsigned long len = theStreamWriter.DefineProperBufferLength();
        std::cout << "\n" << len;
        char* finalBuffer = new char[len];
        memcpy(finalBuffer, &(raw[prevLen]), len);
        std::cout << "\nable to write";
        if (!theStreamWriter.Write(finalBuffer, len)){
            std::cerr << "writing failure:" << "output.dcm" << " at y = " << y << " and z = " << z <<
            std::endl;
            delete [] raw;
            delete [] finalBuffer;
            return 1;
        }
        delete [] finalBuffer;
        prevLen += len;
    }
}
delete raw;

delete[] src; //FIXME

if(dinfo) {
    obj_destroy_decompress(dinfo);
}

obj_image_destroy(image);

return true;
}

bool Different_Resolution( gdcm::StreamImageWriter & theStreamWriter, const char *
    filename, int res, std::ostream& of)
{
    //std::vector<std::string>::const_iterator it = filenames.begin();
    bool b = true;
    int flag = 1;

    gdcm::SmartPointer<gdcm::SequenceOfItems> sq = new
        gdcm::SequenceOfItems();
    sq->SetLengthToUndefined();

    for(int i = res-1 ; i>=0; --i)
    {
        b = b && Write_Resolution( theStreamWriter, filename, i, of ,flag,sq,res);
        // b = b && Get_Resolution( theStreamWriter, filename, i, of ,0);
        flag = 0;
    }
    //b = b && Get_Lowest_Resolution( writer, sq, filename, res-1 );
    //b = b && PopulateSingeFile( writer, sq, jpeg, filename2 );
    //image.SetDimension(2, res )
    return b;
}

```

```

int main(int argc, char *argv[])
{

    if( argc < 4 )
    {
        std::cerr << argv[0] << " input.jp2 output.dcm No. Of Resolutions " << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    char *resolutions = argv[3];
    int res = int((*resolutions)-48);
    //std:: cout << "\nres"<< res;
    gdcm::StreamImageWriter theStreamWriter;

    std::ofstream of;
    of.open( outfile, std::ios::out | std::ios::binary );
    theStreamWriter.SetStream(of);

    if( !Different_Resolution( theStreamWriter, filename,res,of ) ) return 1;

    uint16_t firstTag1 = 0xfffe;
    uint16_t secondTag1 = 0xe0dd;
    uint32_t thirdTag1 = 0x00000000;
    //uint16_t fourthTag1 = 0xffff;
    const int theBufferSize1 = 2*sizeof(uint16_t)+sizeof(uint32_t);
    char* tmpBuffer2 = new char[theBufferSize1];
    memcpy(&(tmpBuffer2[0]), &firstTag1, sizeof(uint16_t));
    memcpy(&(tmpBuffer2[sizeof(uint16_t)]), &secondTag1, sizeof(uint16_t));
    memcpy(&(tmpBuffer2[2*sizeof(uint16_t)]), &thirdTag1, sizeof(uint32_t));
    //memcpy(&(tmpBuffer2[3*sizeof(uint16_t)]), &fourthTag1, sizeof(uint16_t));
    assert( of && !of.eof() && of.good() );
    of.write(tmpBuffer2, theBufferSize1);
    of.flush();
    assert( of );

    return 0;
}

```

## 29.53 ExtractOneFrame.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * This small code shows how to use the gdcm.StreamImageReader API
 * to read a single (whole) frame at a time
 * The API allow extracting a smaller extent of the frame of course.
 * It will write out the extracted frame in /tmp/frame.raw
 *
 * Usage:
 * $ bin/ExtractOneFrame.exe input.dcm
 */
using System;
using gdcm;

public class ExtractOneFrame
{
    public static int Main(string[] args)
    {
        string filename = args[0];

```

```

gdcM.StreamImageReader reader = new gdcM.
    StreamImageReader();

reader.SetFileName( filename );

if (!reader.ReadImageInformation()) return 1;
// Get file infos
gdcM.File f = reader.GetFile();

// get some info about image
UIntArrayType extent = ImageHelper.GetDimensionsValue(f);
//System.Console.WriteLine( extent[0] );
uint dimx = extent[0];
//System.Console.WriteLine( extent[1] );
uint dimy = extent[1];
//System.Console.WriteLine( extent[2] );
uint dimz = extent[2];
PixelFormat pf = ImageHelper.GetPixelFormatValue (f);
int pixelSize = pf.GetPixelSize();
//System.Console.WriteLine( pixelSize );

// buffer to get the pixels
byte[] buffer = new byte[ dimx * dimy * pixelSize ];

for (int i = 0; i < dimz; i++)
{
    // Define that I want the image 0, full size (dimx x dimy pixels)
    reader.DefinePixelExtent(0, (ushort)dimx, 0, (ushort)dimy, (ushort)i, (ushort)(i+1));
    uint buf_len = reader.DefineProperBufferLength(); // take into account pixel size
    //System.Console.WriteLine( buf_len );
    if( buf_len > buffer.Length )
    {
        throw new Exception("buffer is too small for target");
    }

    if (reader.Read(buffer, (uint)buffer.Length))
    {
        using (System.IO.Stream stream =
            System.IO.File.Open(@"tmp/frame.raw",
                System.IO.FileMode.Create))
        {
            System.IO.BinaryWriter writer = new System.IO.BinaryWriter(stream);
            writer.Write(buffer);
        }
    }
    else
    {
        throw new Exception("can't read pixels error");
    }
}

return 0;
}

```

## 29.54 Fake\_Image\_Using\_Stream\_Image\_Writer.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
// This work was realised during the GSOC 2011 by Manoj Alwani

#include "gdcMReader.h"
#include "gdcMMediaStorage.h"
#include "gdcMWriter.h"
#include "gdcMItem.h"

```

```

#include "gdcmImageReader.h"
#include "gdcmAttribute.h"
#include "gdcmFile.h"
#include "gdcmTag.h"
#include "gdcmTransferSyntax.h"
#include "gdcmUIDGenerator.h"
#include "gdcmAnonymizer.h"
#include "gdcmStreamImageWriter.h"
#include "gdcmImageHelper.h"
#include "gdcmTrace.h"

int main(int, char *[])
{
    char * buffer = new char[ 256 * 256 *3 ];
    // *p = (uint8_t*)buffer;
    char * p = buffer;

    gdcm::Trace::DebugOn();
    gdcm::Trace::WarningOn();

    for(int row = 0; row < 256; ++row)
    {
        for(int col = 0; col < 256; ++col)
            //for(int b = 0; b < 256; ++b)
            {
                *p++ = 255;
                *p++ = 0;
                *p++ = 0;
            }

        gdcm::Writer w;
        gdcm::File &file = w.GetFile();
        gdcm::DataSet &ds = file.GetDataSet();

        file.GetHeader().SetDataSetTransferSyntax(
            gdcm::TransferSyntax::ExplicitVRLittleEndian );

        gdcm::UIDGenerator uid;
        gdcm::DataElement de( gdcm::Tag(0x8,0x18) ); // SOP Instance UID
        de.SetVR( gdcm::VR::UI );
        const char *u = uid.Generate();
        de.SetByteValue( u, strlen(u) );
        ds.Insert( de );

        gdcm::DataElement del( gdcm::Tag(0x8,0x16) );
        del.SetVR( gdcm::VR::UI );
        gdcm::MediaStorage ms(
            gdcm::MediaStorage::VLWholeSlideMicroscopyImageStorage
        );
        del.SetByteValue( ms.GetString(), strlen(ms.GetString()) );
        ds.Insert( del );

        const char mystr[] = "RGB";
        gdcm::DataElement de2( gdcm::Tag(0x28,0x04) );
        //de.SetTag(gdcm::Tag(0x28,0x04));
        de2.SetVR( gdcm::VR::CS );
        de2.SetByteValue(mystr, strlen(mystr));
        ds.Insert( de2 );

        gdcm::Attribute<0x0028,0x0010> row = {256};
        //row.SetValue(512);
        ds.Insert( row.GetAsDataElement() );
        // w.SetCheckFileMetaInformation( true );
        gdcm::Attribute<0x0028,0x0011> col = {256};
        ds.Insert( col.GetAsDataElement() );

        gdcm::Attribute<0x0028,0x0008> Number_Of_Frames = {1};
        ds.Insert( Number_Of_Frames.GetAsDataElement() );

        gdcm::Attribute<0x0028,0x0100> at = {8};
        ds.Insert( at.GetAsDataElement() );

        gdcm::Attribute<0x0028,0x0002> at1 = {3}; //bits per pixel
        ds.Insert( at1.GetAsDataElement() );

        gdcm::Attribute<0x0028,0x0101> at2 = {8};
        ds.Insert( at2.GetAsDataElement() );

        gdcm::Attribute<0x0028,0x0102> at3 = {7};
    }
}

```

```

    ds.Insert( at3.GetAsDataElement() );

    gdcmm::Attribute<0x0028,0x006> at4 = {0};
    ds.Insert( at4.GetAsDataElement() );

    gdcmm::Attribute<0x0028,0x0103> at5 = {0};
    ds.Insert( at5.GetAsDataElement() );

    //de.SetTag(gdcmm::Tag(0x7fe0,0x0010));
    //ds.Insert(de);

    gdcmm::StreamImageWriter theStreamWriter;
    gdcmm::SmartPointer<gdcmm::SequenceOfItems> sq = new
        gdcmm::SequenceOfItems();
    sq->SetLengthToUndefined();

    uint16_t row1 = 256;
    uint16_t col1 = 256;
    //std::cout << row;

    gdcmm::Element<gdcmm::VR::IS,gdcmm::VM::VM1> el2;
    el2.SetValue(1);
    gdcmm::DataElement rfn = el2.GetAsDataElement();    //rfn --->
        reference frame number
    rfn.SetTag( gdcmm::Tag(0x0008,0x1160) );

    gdcmm::Element<gdcmm::VR::US,gdcmm::VM::VM2> el;
    el.SetValue(1,0);
    el.SetValue(1,1);
    gdcmm::DataElement ulr = el.GetAsDataElement();    //ulr --> upper
        left col/row
    ulr.SetTag( gdcmm::Tag(0x0048,0x0201) );

    gdcmm::Element<gdcmm::VR::US,gdcmm::VM::VM2> el1;
    el1.SetValue(col1,0);
    el1.SetValue(row1,1);
    gdcmm::DataElement brr = el1.GetAsDataElement();
    brr.SetTag( gdcmm::Tag(0x0048,0x0202) );    //brr --> bottom right col/row

    gdcmm::Item it;
    gdcmm::DataSet &nds = it.GetNestedDataSet();
    nds.Insert( rfn );
    nds.Insert(ulr);
    nds.Insert(brr);

    sq->AddItem(it);

    gdcmm::DataElement des( gdcmm::Tag(0x0048,0x0200) );
    des.SetVR(gdcmm::VR::SQ);
    des.SetValue(*sq);
    des.SetVLToUndefined();

    ds.Insert( des );

    theStreamWriter.SetFile(file);

    std::ofstream of;
    of.open( "output.dcm", std::ios::out | std::ios::binary );
    theStreamWriter.SetStream(of);

    if (!theStreamWriter.CanWriteFile()){
        delete [] buffer;
        std::cout << "Not able to write";
        return 0; //this means that the file was unwritable, period.
        //very similar to a ReadImageInformation failure
    }
    else
        std::cout<<"\nable to read";

    if (!theStreamWriter.WriteImageInformation()){
        std::cerr << "unable to write image information" << std::endl;
        delete [] buffer;
        return 1; //the CanWrite function should prevent getting here, else,
        //that's a test failure
    }

    std::vector<unsigned int> extent =
        gdcmm::ImageHelper::GetDimensionsValue(file);

```

```

unsigned short xmax = extent[0];
unsigned short ymax = extent[1];
unsigned short theChunkSize = 1;
unsigned short ychunk = extent[1]/theChunkSize; //go in chunk sizes of theChunkSize
unsigned short zmax = extent[2];

std::cout << xmax << ymax << zmax;

if (xmax == 0 || ymax == 0)
{
    std::cerr << "Image has no size, unable to write zero-sized image." << std::endl;
    return 0;
}

int z, y, nexty;
unsigned long prevLen = 0; //when going through the char buffer, make sure to grab
//the bytes sequentially. So, store how far you got in the buffer with each iteration.
for (z = 0; z < zmax; ++z){
    for (y = 0; y < ymax; y += ychunk){
        nexty = y + ychunk;
        if (nexty > ymax) nexty = ymax;
        theStreamWriter.DefinePixelExtent(0, xmax, y, nexty, z, z+1);
        unsigned long len = theStreamWriter.DefineProperBufferLength();
        std::cout << "\n" << len;
        char* finalBuffer = new char[len];
        memcpy(finalBuffer, &(buffer[prevLen]), len);
        std::cout << "\nable to write";
        if (!theStreamWriter.Write(finalBuffer, len)){
            std::cerr << "writing failure:" << "output.dcm" << " at y = " << y << " and z = " << z <<
            std::endl;
            delete [] buffer;
            delete [] finalBuffer;
            return 1;
        }
        delete [] finalBuffer;
        prevLen += len;
    }
}
delete buffer;

uint16_t firstTag1 = 0xfffe;
uint16_t secondTag1 = 0xe0dd;
uint32_t thirdTag1 = 0x00000000;
//uint16_t fourthTag1 = 0xffff;
const int theBufferSize1 = 2*sizeof(uint16_t)+sizeof(uint32_t);
char* tmpBuffer2 = new char[theBufferSize1];
memcpy(&(tmpBuffer2[0]), &firstTag1, sizeof(uint16_t));
memcpy(&(tmpBuffer2[sizeof(uint16_t)]), &secondTag1, sizeof(uint16_t));
memcpy(&(tmpBuffer2[2*sizeof(uint16_t)]), &thirdTag1, sizeof(uint32_t));
//memcpy(&(tmpBuffer2[3*sizeof(uint16_t)]), &fourthTag1, sizeof(uint16_t));
assert( of && !of.eof() && of.good() );
of.write(tmpBuffer2, theBufferSize1);
of.flush();
assert( of );

return 0;
}

```

## 29.55 FileAnonymize.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
* Simple C# example

```

```

*
* Usage:
* $ mono bin/FileAnonymize.exe input.dcm output.dcm
*/
using System;
using gdcm;

public class FileAnonymize
{
    public static int Main(string[] args)
    {
        string filename = args[0];
        string outfilename = args[1];

        gdcm.FileAnonymizer fa = new gdcm.FileAnonymizer();
        fa.SetInputFileName( filename );
        fa.SetOutputFileName( outfilename );

        // Empty Operations
        // It will create elements, since those tags are non-registered public elements (2011):
        fa.Empty( new Tag(0x0008,0x1313) );
        fa.Empty( new Tag(0x0008,0x1317) );
        // Remove Operations
        // The following Tag are actually carefully chosen, since they refer to SQ:
        fa.Remove( new Tag(0x0008,0x2112) );
        fa.Remove( new Tag(0x0008,0x9215) );
        // Replace Operations
        // do not call replace operation on SQ attribute !
        fa.Replace( new Tag(0x0018,0x5100), "MYVALUE " );
        fa.Replace( new Tag(0x0008,0x1160), "MYOTHERVAL" );

        if ( !fa.Write() )
        {
            System.Console.WriteLine( "Could not write" );
            return 1;
        }

        return 0;
    }
}

```

## 29.56 FileAnonymize.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

import gdcm.*;

public class FileAnonymize
{
    public static class MyWatcher extends SimpleSubjectWatcher
    {
        public MyWatcher(Subject s) { super(s,"Override String"); }
        protected void ShowProgress(Subject caller, Event evt)
        {
            ProgressEvent pe = ProgressEvent.Cast(evt);
            System.out.println( "This is my progress: " + pe.GetProgress() );
        }
    }

    public static void main(String[] args) throws Exception
    {
        String input = args[0];
        String output = args[1];

        FileAnonymizer fa = new FileAnonymizer();
    }
}

```

```

fa.SetInputFileName( input );
fa.SetOutputFileName( output );

// Empty Operations
// It will create elements, since those tags are non-registered public elements (2011):
fa.Empty( new Tag(0x0008,0x1313) );
fa.Empty( new Tag(0x0008,0x1317) );
// Remove Operations
// The following Tag are actually carefully chosen, since they refer to SQ:
fa.Remove( new Tag(0x0008,0x2112) );
fa.Remove( new Tag(0x0008,0x9215) );
// Replace Operations
// do not call replace operation on SQ attribute !
fa.Replace( new Tag(0x0018,0x5100), "MYVALUE " );
fa.Replace( new Tag(0x0008,0x1160), "MYOTHERVAL" );

if( !fa.Write() )
{
    System.out.println( "Could not write" );
    return;
}

System.out.println( "success" );
}
}

```

## 29.57 FileChangeTS.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Simple C# example
 *
 * Shows multiple steps:
 * Steps 1.
 * Create a fake (dummy) DICOM file, with size 512 x 512 x 2 We use a small
 * image to be able to create the volume in memory Of course you can use any
 * existing DICOM instead
 *
 * Step 2.
 * Hack the DICOM file to pretend the number of frames is 1000 (instead of 2)
 * At this point in time this makes the DICOM file invalid (truncated). But the
 * next step will fix this.
 *
 * Step 3.
 * Use C# to create a binary data which will represent our source object for
 * image.
 *
 * Step 4.
 * We use gdcm.FileStreamer to merge the template DICOM file from Step 2, with
 * the binary data from Step 3. We decide to read a scanline at a time, but
 * this can be read with any number of bytes. AppendToDataElement() will always
 * do the proper computation.
 *
 * Step 5.
 * We compress this gigantic file, into [JPEG Lossless, Non-Hierarchical,
 * First-Order Prediction (Process 14 [Selection Value 1])]
 *
 * Usage:
 * $ mono bin/FileChangeTS.exe small.dcm big.dcm raw.data merge.dcm jpeg.dcm
 */
using System;
using System.IO;
using gdcm;

```



```

public class FileChangeTS
{
    public static byte[] StrToByteArray(string str)
    {
        System.Text.ASCIIEncoding encoding=new System.Text.ASCIIEncoding();
        return encoding.GetBytes(str);
    }
    // Create a 256 x 256 Secondary Capture Image Storage
    static private void CreateSmallDICOM(string fileName)
    {
        using( var writer = new gdcm.PixmapWriter() )
        {
            gdcm.Pixmap img = writer.GetImage();
            img.SetNumberOfDimensions( 3 );
            img.SetDimension(0, 512 );
            img.SetDimension(1, 512 );
            img.SetDimension(2, 2 ); // fake a 3d volume
            PhotometricInterpretation pi = new PhotometricInterpretation( PhotometricInterpretation.PIType.
            MONOCHROME2 );
            img.SetPhotometricInterpretation( pi );
            gdcm.DataElement pixeldata = new gdcm.DataElement( new
            gdcm.Tag(0x7fe0,0x0010) );
            byte[] buffer = new byte[ 512 * 512 * 2 ];
            pixeldata.SetByteValue( buffer, new gdcm.VL((uint)buffer.Length) );
            img.SetDataElement( pixeldata );

            gdcm.File file = writer.GetFile();
            gdcm.DataSet ds = file.GetDataSet();
            gdcm.DataElement ms = new gdcm.DataElement(new
            gdcm.Tag(0x0008,0x0016));
            string mediastorage = "1.2.840.10008.5.1.4.1.1.7.2"; // Multi-frame Grayscale Byte Secondary Capture
            Image Storage
            byte[] val = StrToByteArray(mediastorage);
            ms.SetByteValue( val, new gdcm.VL( (uint)val.Length) );
            ds.Insert( ms );

            writer.SetFileName( fileName );
            writer.Write();
        }
    }
    static private void CreateBigDICOM(string fileName, string outfilename)
    {
        using( var ano = new gdcm.FileAnonymizer() )
        {
            // The following is somewhat dangerous, do not try at home:
            string nframes = "1000";
            ano.Replace( new gdcm.Tag(0x0028,0x0008), nframes );
            ano.SetInputFileName(fileName);
            ano.SetOutputFileName(outfilename);
            ano.Write(); // at this point the DICOM is invalid !
        }
    }
    static private void CreateDummyFile(string fileName, long length)
    {
        using (var fileStream = new FileStream(fileName, FileMode.Create, FileAccess.Write, FileShare.None))
        {
            // Looks like C# always init to 0 (fallocate ?)
            // For the purpose of the test we could add some random noise
            fileStream.SetLength(length);
        }
    }
    static private void ReadBytesIntoArray( byte[] array, FileStream source )
    {
        int numBytesToRead = array.Length;
        int numBytesRead = 0;
        while (numBytesToRead > 0)
        {
            // According to spec: Read() may return anything from 0 to numBytesToRead.
            int n = source.Read(array, numBytesRead, numBytesToRead);

            // Break when the end of the file is reached.
            if (n == 0)
                break;

            numBytesRead += n;
            numBytesToRead -= n;
        }
    }
    static private void AssembleDICOMAndRaw(string dicomfn, string rawdata, string outfn)
    {
        using ( var fs = new gdcm.FileStreamer() )
    }

```

```

    {
        fs.SetTemplateFileName(dicomfn);
        fs.SetOutputFileName(outfn);
        gdcm.Tag pixeldata = new gdcm.Tag(0x7fe0, 0x0010);
        // FileStreamer support automatic checking of pixel data length
        // based on DICOM attributes, only if we say so:
        fs.CheckDataElement( pixeldata );
        // Declare we are working on Pixel Data attribute:
        fs.StartDataElement( pixeldata );
        using (FileStream rawSource = new FileStream(rawdata,
            FileMode.Open, FileAccess.Read))
        {
            byte[] bytes = new byte[512];
            // Only read one scanline at a time
            // We could have been reading more at once, if this is more efficient,
            // AppendToDataElement will do the logic in all cases.
            for( int i = 0; i < 512 * 1000; ++i )
            {
                // Read the source file into a byte array.
                ReadBytesIntoArray( bytes, rawSource );
                fs.AppendToDataElement( pixeldata, bytes, (uint)bytes.Length );
            }
        }
        if( !fs.StopDataElement( pixeldata ) )
        {
            // Most likely an issue with Pixel Data Length computation:
            throw new Exception("StopDataElement failed");
        }
    }
}

static private void CompressIntoJPEG(string rawdicom, string jpegdicom)
{
    using( var sfcts = FileChangeTransferSyntax.New() )
    {
        // Need to retrieve the actual C++ reference, to pass to
        // SimpleSubjectWatcher:
        FileChangeTransferSyntax fcts = sfcts.__ref__();
        SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(fcts, "FileChangeTransferSyntax");
        gdcm.TransferSyntax ts = new TransferSyntax( TransferSyntax.TSType.
            JPEGLosslessProcess14_1 );
        fcts.SetTransferSyntax( ts );
        fcts.SetInputFileName( rawdicom );
        fcts.SetOutputFileName( jpegdicom );
        fcts.Change();
    }
}

public static int Main(string[] args)
{
    string filename = args[0];
    string outfilename = args[1];
    string rawfilename = args[2];
    string mergefn = args[3];
    string jpegfn = args[4];

    CreateSmallDICOM(filename);
    CreateBigDICOM(filename, outfilename);
    CreateDummyFile(rawfilename, 512 * 512 * 1000 );
    AssembleDICOMAndRaw(outfilename, rawfilename, mergefn);
    CompressIntoJPEG(mergefn, jpegfn);

    return 0;
}
}

```

## 29.58 FileChangeTSLossy.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====

```

```

=====*/

/*
 * Simple C# example
 *
 * Shows multiple steps:
 * Steps 1.
 * Create a fake (dummy) DICOM file, with size 512 x 512 x 2 We use a small
 * image to be able to create the volume in memory Of course you can use any
 * existing DICOM instead
 *
 * Step 2.
 * Hack the DICOM file to pretend the number of frames is 1000 (instead of 2)
 * At this point in time this makes the DICOM file invalid (truncated). But the
 * next step will fix this.
 *
 * Step 3.
 * Use C# to create a binary data which will represent our source object for
 * image.
 *
 * Step 4.
 * We use gdcm.FileStreamer to merge the template DICOM file from Step 2, with
 * the binary data from Step 3. We decide to read a scanline at a time, but
 * this can be read with any number of bytes. AppendToDataElement() will always
 * do the proper computation.
 *
 * Step 5.
 * We compress this gigantic file, into [JPEG Baseline (Process 1): Default Transfer Syntax for Lossy JPEG
 * 8 Bit Image Compression]
 *
 * Usage:
 * $ bin/FileChangeTSLossy.exe small.dcm big.dcm raw.data merge.dcm jpeg.dcm
 */
using System;
using System.IO;
using gdcm;

public class FileChangeTS
{
    public static byte[] StrToByteArray(string str)
    {
        System.Text.ASCIIEncoding encoding=new System.Text.ASCIIEncoding();
        return encoding.GetBytes(str);
    }
    // Create a 256 x 256 Secondary Capture Image Storage
    static private void CreateSmallDICOM(string fileName)
    {
        using( var writer = new gdcm.PixmapWriter() )
        {
            gdcm.Pixmap img = writer.GetImage();
            img.SetNumberOfDimensions( 3 );
            img.SetDimension(0, 512 );
            img.SetDimension(1, 512 );
            img.SetDimension(2, 2 ); // fake a 3d volume
            PhotometricInterpretation pi = new PhotometricInterpretation( PhotometricInterpretation.PIType.
MONOCHROME2 );
            img.SetPhotometricInterpretation( pi );
            gdcm.DataElement pixeldata = new gdcm.DataElement( new
gdcm.Tag(0x7fe0,0x0010) );
            byte[] buffer = new byte[ 512 * 512 * 2 ];
            pixeldata.SetByteValue( buffer, new gdcm.VL((uint)buffer.Length) );
            img.SetDataElement( pixeldata );

            gdcm.File file = writer.GetFile();
            gdcm.DataSet ds = file.GetDataSet();
            gdcm.DataElement ms = new gdcm.DataElement(new
gdcm.Tag(0x0008,0x0016));
            string mediastorage = "1.2.840.10008.5.1.4.1.1.7.2"; // Multi-frame Grayscale Byte Secondary Capture
Image Storage
            byte[] val = StrToByteArray(mediastorage);
            ms.SetByteValue( val, new gdcm.VL( (uint)val.Length) );
            ds.Insert( ms );

            writer.SetFileName( fileName );
            writer.Write();
        }
    }
    static private void CreateBigDICOM(string fileName, string outfilename)
    {
        using( var ano = new gdcm.FileAnonymizer() )

```

```

    {
        // The following is somewhat dangerous, do not try at home:
        string nframes = "1000";
        ano.Replace( new gdcm.Tag(0x0028,0x0008), nframes );
        ano.SetInputFileName(fileName);
        ano.SetOutputFileName(outfilename);
        ano.Write(); // at this point the DICOM is invalid !
    }
}

static private void CreateDummyFile(string fileName, long length)
{
    using (var fileStream = new FileStream(fileName, FileMode.Create, FileAccess.Write, FileShare.None))
    {
        // Looks like C# always init to 0 (fallocate ?)
        // For the purpose of the test we could add some random noise
        fileStream.SetLength(length);
    }
}

static private void ReadBytesIntoArray( byte[] array, FileStream source )
{
    int numBytesToRead = array.Length;
    int numBytesRead = 0;
    while (numBytesToRead > 0)
    {
        // According to spec: Read() may return anything from 0 to numBytesToRead.
        int n = source.Read(array, numBytesRead, numBytesToRead);

        // Break when the end of the file is reached.
        if (n == 0)
            break;

        numBytesRead += n;
        numBytesToRead -= n;
    }
}

static private void AssembleDICOMAndRaw(string dicomfn, string rawdata, string outfn)
{
    using ( var fs = new gdcm.FileStreamer() )
    {
        fs.SetTemplateFileName(dicomfn);
        fs.SetOutputFileName(outfn);
        gdcm.Tag pixeldata = new gdcm.Tag(0x7fe0, 0x0010);
        // FileStreamer support automatic checking of pixel data length
        // based on DICOM attributes, only if we say so:
        fs.CheckDataElement( pixeldata );
        // Declare we are working on Pixel Data attribute:
        fs.StartDataElement( pixeldata );
        using (FileStream rawSource = new FileStream(rawdata,
            FileMode.Open, FileAccess.Read))
        {
            byte[] bytes = new byte[512];
            // Only read one scanline at a time
            // We could have been reading more at once, if this is more efficient,
            // AppendToDataElement will do the logic in all cases.
            for( int i = 0; i < 512 * 1000; ++i )
            {
                // Read the source file into a byte array.
                ReadBytesIntoArray( bytes, rawSource );
                fs.AppendToDataElement( pixeldata, bytes, (uint)bytes.Length );
            }
        }
        if( !fs.StopDataElement( pixeldata ) )
        {
            // Most likely an issue with Pixel Data Length computation:
            throw new Exception("StopDataElement failed");
        }
    }
}

static private void CompressIntoJPEG(string rawdicom, string jpegdicom)
{
    using( var sfcts = FileChangeTransferSyntax.New() )
    {
        // Need to retrieve the actual C++ reference, to pass to
        // SimpleSubjectWatcher:
        FileChangeTransferSyntax fcts = sfcts.__ref__();
        SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(fcts, "FileChangeTransferSyntax");
        gdcm.TransferSyntax ts = new TransferSyntax( TransferSyntax.TSType.
            JPEGBaselineProcess1 );
        fcts.SetTransferSyntax( ts );
        ImageCodec ic = fcts.GetCodec();
        JPEGCodec jpeg = JPEGCodec.Cast( ic );
    }
}

```

```

        jpeg.SetLossless( false );
        jpeg.SetQuality( 50 ); // poor quality !

        fcts.SetInputFileName( rawdicom );
        fcts.SetOutputFileName( jpegdicom );
        fcts.Change();
    }
}

public static int Main(string[] args)
{
    string filename = args[0];
    string outfilename = args[1];
    string rawfilename = args[2];
    string mergefn = args[3];
    string jpegfn = args[4];

    CreateSmallDICOM(filename);
    CreateBigDICOM(filename, outfilename);
    CreateDummyFile(rawfilename, 512 * 512 * 1000 );
    AssembleDICOMAndRaw(outfilename, rawfilename, mergefn);
    CompressIntoJPEG(mergefn, jpegfn);

    return 0;
}
}

```

## 29.59 FileStreaming.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Simple C# example
 *
 * Usage:
 * $ mono bin/FileStreaming.exe gdcmData/CT_16b_signed-UsedBits13.dcm output.dcm
 *
 * The class will take care of group handling and will use the first available group:
 * (0009,0012) ?? (LO) [MYTEST] # 6,1 Private Creator
 */
using System;
using gdcm;

public class FileStreaming
{
    public static int Main(string[] args)
    {
        string filename = args[0];
        string outfilename = args[1];

        gdcm.PrivateTag pt = new gdcm.PrivateTag( new
            gdcm.Tag(0x9,0x10), "MYTEST" );

        gdcm.FileStreamer fs = new gdcm.FileStreamer();
        fs.SetTemplateFileName( filename );
        fs.SetOutputFileName( outfilename );

        byte[] buffer = new byte[ 8192 ];
        uint len = (uint)buffer.Length;

        // In this example, we want that each newly created Private Attribute
        // contains at most 1000 bytes of incoming dataset.
        // We are also calling the function twice to check that appending mode is
        // working from one call to the other. The last element will have a length
        // of (2 * 8192) % 1000 = 384
        if( !fs.StartGroupDataElement( pt, 1000, 1 )

```

```

    || !fs.AppendToGroupDataElement( pt, buffer, len )
    || !fs.AppendToGroupDataElement( pt, buffer, len )
    || !fs.StopGroupDataElement( pt ) )
    {
        System.Console.WriteLine( "Could not change private group" );
        return 1;
    }

    return 0;
}
}

```

## 29.60 FindAllPatientName.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14 """
15 This example shows how one can use the gdcm.CompositeNetworkFunctions class
16 for executing a C-FIND query
17 It will print the list of patient name found
18
19 Usage:
20
21 python FindAllPatientName.py
22
23 """
24
25 import gdcm
26
27 # Patient Name
28 tag = gdcm.Tag(0x10,0x10)
29 de = gdcm.DataElement(tag)
30
31 # Search all patient name where string match 'F*'
32 de.SetByteValue('F*',gdcm.VL(2))
33
34 ds = gdcm.DataSet()
35 ds.Insert(de)
36
37 cnf = gdcm.CompositeNetworkFunctions()
38 theQuery = cnf.ConstructQuery(gdcm.ePatientRootType,gdcm.ePatient,ds)
39
40 #print theQuery.ValidateQuery()
41
42 # prepare the variable for output
43 ret = gdcm.DataSetArrayType()
44
45 # Execute the C-FIND query
46 cnf.CFind('dicom.example.com',11112,theQuery,ret,'GDCM_PYTHON','ANY-SCP')
47
48 for i in range(0,ret.size()):
49     print "Patient #",i
50     print ret[i]

```

## 29.61 FixBrokenJ2K.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.

```

See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the above copyright notice for more information.

```

=====*/
#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmImageReader.h"
#include "gdcmSequenceOfFragments.h"
#include "gdcmFile.h"

// http://www.lost.in.ua/dicom/c.dcm
//
// -> BuggyJ2Kvvvua-fixed2-j2k.dcm

/*
 * This program attempts to fix a broken J2K/DICOM:
 * It contains 2 bugs:
 * 1. The first 8 bytes seems to be random bytes: remove them
 * 2. YCC is set to 1, while image is grayscale need to set it back to 0
 *
 * Ref:
 * It's a software from http://rentgenprom.ru/ , shipped with universal digital radiographic units
 * "ProScan-2000". The Ukrainian manufacturer developed own digital radiographic unit and it is
 * compatible with software from "ProScan-2000".
 */
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::File &file = reader.GetFile();
    const gdcm::DataElement &pixeldata0 = file.GetDataSet().
        GetDataElement( gdcm::Tag(0x7fe0,0x0010) );
    const gdcm::SequenceOfFragments *sqf = pixeldata0.
        GetSequenceOfFragments();
    if( !sqf )
    {
        return 1;
    }
    const gdcm::Fragment &frag0 = sqf->GetFragment(0);

    const gdcm::ByteValue *bv = frag0.GetByteValue();
    const char *ptr = bv->GetPointer();
    size_t len = bv->GetLength();

    const char sig[] = "\x00\x00\x00\x00\x6A\x70\x32\x63";
    if( memcmp(ptr, sig, sizeof(sig)) != 0 )
    {
        std::cerr << "magic random signature not found" << std::endl;
        return 1;
    }

    // Apparently the flag to enable a color transform on 3 color components is set in
    // the COD marker. (YCC is byte[6] in the COD marker)
    // we need to disable this flag;
    const char *cod_marker = ptr + 0x35; /* 0x2d + 0x8 */ // FIXME
    if( cod_marker[0] == (char)0xff && cod_marker[1] == 0x52 )
    {
        // found start of COD
        if( cod_marker[6+2] == 1 )
        {
            // Change in place:
            *((char*)cod_marker + 6+2) = 0;
            // Prepare a new DataElement:
            gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
            pixeldata.SetVR( gdcm::VR::OB );
            gdcm::SmartPointer<gdcm::SequenceOfFragments> sq = new

```

```

        gdcm::SequenceOfFragments;

        gdcm::Fragment frag;
        // remove 8 first bytes:
        frag.SetByteValue( ptr + 8, (uint32_t)(len - 8) );
        sq->AddFragment( frag );
        pixeldata.SetValue( *sq );
        file.GetDataSet().Replace( pixeldata );
    }
    else
    {
        return 1;
    }
}
else
{
    std::cerr << "COD not found" << (int)cod_marker[0] << std::endl;
    return 1;
}

gdcm::Writer writer;
writer.SetFile( reader.GetFile() );
writer.SetFileName( outfilename );
writer.CheckFileMetaInformationOff();
if( !writer.Write() )
{
    std::cerr << "Could not write" << std::endl;
}

// paranoid check:
gdcm::ImageReader ireader;
ireader.SetFileName( outfilename );
if( !ireader.Read() )
{
    std::cerr << "file written is still not valid, please report" << std::endl;
    return 1;
}

return 0;
}

```

## 29.62 FixCommaBug.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Using LC_NUMERIC set to something not compatible with "C" it is possible to write out "," instead of
17 "." as required by the DICOM standard
18 Issue is still current (IMHO) with gdcm 2.0.9
19 """
20
21 import gdcm
22 import sys
23
24 filename = sys.argv[1]
25 outname = sys.argv[2]
26
27 # read
28 r = gdcm.Reader()
29 r.SetFileName( filename )
30 if not r.Read():
31     print "not valid"
32     sys.exit(1)
33

```



```

34 file = r.GetFile()
35 dataset = file.GetDataSet()
36
37 ano = gdcm.Anonymizer()
38 ano.SetFile( file )
39
40 tags = [
41     gdcm.Tag(0x0018,0x1164),
42     gdcm.Tag(0x0018,0x0088),
43     gdcm.Tag(0x0018,0x0050),
44     gdcm.Tag(0x0028,0x0030),
45 ]
46
47 for tag in tags:
48     print tag
49     if dataset.FindDataElement( tag ):
50         pixelpacing = dataset.GetDataElement( tag )
51         #print pixelpacing
52         bv = pixelpacing.GetByteValue()
53         str = bv.GetBuffer()
54         #print bv.GetLength()
55         #print len(str)
56         new_str = str.replace(",",".")
57         # Need to explicitly pass bv.GetLength() to remove any trailing garbage
58         ano.Replace( tag, new_str, bv.GetLength() )
59
60 #print dataset
61
62 w = gdcm.Writer()
63 w.SetFile( file )
64 w.SetFileName( outname )
65 if not w.Write():
66     print "Cannot write"
67     sys.exit(1)
68
69 # paranoid:
70 image_reader = gdcm.ImageReader()
71 image_reader.SetFileName( outname )
72 if not image_reader.Read():
73     print "there is still a comma"
74     sys.exit(1)
75
76 print "Success!"
77 sys.exit(0) # success

```

## 29.63 FixJAIBugJPEGLS.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmImageReader.h"

#include <fstream>

#include "gdcm_charls.h"

/*
 * This small example should show how one can handle the famous JAI-JPEGLS bug
 * It will take in as invalid DICOM/JAI-JPEG-LS and write out as Explicit Little
 * Endian. One can use 'gdcmconv --jpegls' to recompress properly
 *
 * References:
 * http://charls.codeplex.com/discussions/230307?ProjectName=charls
 * http://charls.codeplex.com/workitem/7297
 * http://www.dcm4che.org/jira/browse/DCM-442

```

```

* http://www.dcm4che.org/jira/browse/DCMEE-1144
* http://java.net/jira/browse/JAI_IMAGEIO_CORE-183
*
* Explanation of the issue:
*
* Seems, the error is in the calculation of the default values for thresholds T1,
* T2, T3, in particular min(MAXVAL, 4095) is not applied in
*
* FACTOR = (min(MAXVAL, 4095) + 128)/256
*
* as specified in http://www.itu.int/rec/T-REC-T.87-199806-I/en .
*
*/
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::FileMetaInformation::SetSourceApplicationEntityTitle
        ( "FixJAIBugJPEGs" );

    gdcm::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::Image &image = reader.GetImage();
    //unsigned long len = image.GetBufferLength();
    const gdcm::DataElement &in =
        reader.GetFile().GetDataSet().GetDataElement(
            gdcm::Tag(0x7fe0,0x0010) );
    const gdcm::SequenceOfFragments *sf = in.
        GetSequenceOfFragments();
    if( !sf )
    {
        std::cerr << "No pixel data (or not encapsulated)" << std::endl;
        return 1;
    }
    const unsigned int *dims = image.GetDimensions();
    if ( sf->GetNumberOfFragments() != dims[2] )
    {
        std::cerr << "Unsupported" << std::endl;
        return 1;
    }

    // unsigned long totalLen = sf->ComputeByteLength();
    std::vector<BYTE> rbyteOutall;
    for(unsigned int i = 0; i < sf->GetNumberOfFragments(); ++i)
    {
        const gdcm::Fragment &frag = sf->GetFragment(i);
        if( frag.IsEmpty() ) return 1;
        const gdcm::ByteValue *bv = frag.GetByteValue();
        if( !bv ) return 1;
        unsigned long totalLen = bv->GetLength();

        std::vector<char> vbuffer;
        vbuffer.resize( totalLen );
        char *buffer = &vbuffer[0];
        bv->GetBuffer(buffer, totalLen);
        const BYTE* pbyteCompressed0 = (const BYTE*)buffer;
        while( totalLen > 0 && pbyteCompressed0[totalLen-1] != 0xd9 )
        {
            totalLen--;
        }

        JlsParameters metadata;
        if (JpegLsReadHeader(buffer, totalLen, &metadata) != OK)
        {
            std::cerr << "Cant parse jpegls" << std::endl;
            return false;
        }

        std::cout << metadata.width << std::endl;
        std::cout << metadata.height << std::endl;
        std::cout << metadata.bitspersample << std::endl;
    }
}

```

```

gdcM::PixelFormat const & pf = image.GetPixelFormat();
std::cout << pf << std::endl;

// http://charls.codeplex.com/discussions/230307?ProjectName=charls
unsigned char marker_lse_13[] = {
    0xFF, 0xF8, 0x00, 0x0D,
    0x01,
    0x1F, 0xFF,
    0x00, 0x22, // T1 = 34
    0x00, 0x83, // T2 = 131
    0x02, 0x24, // T3 = 548
    0x00, 0x40
};

unsigned char marker_lse_14[] = {
    0xFF, 0xF8, 0x00, 0x0D,
    0x01,
    0x3F, 0xFF,
    0x00, 0x42, // T1 = 66
    0x01, 0x03, // T2 = 259
    0x04, 0x44, // T3 = 1092
    0x00, 0x40
};

unsigned char marker_lse_15[] = {
    0xFF, 0xF8, 0x00, 0x0D,
    0x01,
    0x7F, 0xFF,
    0x00, 0x82, // T1 = 130
    0x02, 0x03, // T2 = 515
    0x08, 0x84, // T3 = 2180
    0x00, 0x40
};

unsigned char marker_lse_16[] = {
    0xFF, 0xF8, 0x00, 0x0D,
    0x01,
    0xFF, 0xFF,
    0x01, 0x02, // T1 = 258
    0x04, 0x03, // T2 = 1027
    0x11, 0x04, // T3 = 4356
    0x00, 0x40
};

const unsigned char *marker_lse = NULL;
switch( metadata.bitspersample )
{
case 13:
    marker_lse = marker_lse_13;
    break;
case 14:
    marker_lse = marker_lse_14;
    break;
case 15:
    marker_lse = marker_lse_15;
    break;
case 16:
    marker_lse = marker_lse_16;
    break;
}
if( !marker_lse )
{
    std::cerr << "Cant handle: " << metadata.bitspersample << std::endl;
    return 1;
}

// FIXME: One should recompute the value for 0x0F
vbuffer.insert( vbuffer.begin() + 0x0F, marker_lse, marker_lse+15);

#if 0
std::ofstream of( "/tmp/d.jls", std::ios::binary );
of.write( &vbuffer[0], vbuffer.size() );
of.close();
#endif

const char *pbyteCompressed = &vbuffer[0];
size_t cbyteCompressed = vbuffer.size(); // updated legnth

JlsParameters params;
JpegLsReadHeader(pbyteCompressed, cbyteCompressed, &params);

```

```

std::vector<BYTE> rgbyteOut;
//rgbyteOut.resize( image.GetBufferLength() );
rgbyteOut.resize(params.height *params.width * ((params.bitspersample + 7)
/ 8) * params.components);

JLS_ERROR result =
JpegLsDecode(&rgbyteOut[0], rgbyteOut.size(), pbyteCompressed, cbyteCompressed, &params );
if (result != OK)
{
std::cerr << "Could not patch JAI-JPEGLS" << std::endl;
return 1;
}
rgbyteOutall.insert( rgbyteOutall.end(), rgbyteOut.begin(), rgbyteOut.end() );
}

gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
pixeldata.SetVR( gdcm::VR::OW );
pixeldata.SetByteValue( (char*)&rgbyteOutall[0], (uint32_t)rgbyteOutall.size() );

// Add the pixel data element
reader.GetFile().GetDataSet().Replace( pixeldata );
reader.GetFile().GetHeader().SetDataSetTransferSyntax(
gdcm::TransferSyntax::ExplicitVRLittleEndian);

gdcm::Writer writer;
writer.SetFileName( outfilename );
writer.SetFile( reader.GetFile() );
writer.Write();

std::cout << "Success !" << std::endl;

return 0;
}

```

## 29.64 gdcmorthoplanes.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

#include "vtkActor.h"
#include "vtkCamera.h"
#include "vtkMatrix4x4.h"
#include "vtkTransform.h"
#include "vtkAssembly.h"
#include "vtkCellPicker.h"
#include "vtkCommand.h"
#include "vtkImageActor.h"
#include "vtkImageMapToColors.h"
#include "vtkImageOrthoPlanes.h"
#include "vtkImagePlaneWidget.h"
#include "vtkImageReader.h"
#include "vtkInteractorEventRecorder.h"
#include "vtkLookupTable.h"
#include "vtkOutlineFilter.h"
#include "vtkPolyDataMapper.h"
#include "vtkProperty.h"
#include "vtkRenderWindow.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkRenderer.h"
#include "vtkVolume16Reader.h"
#include "vtkImageData.h"
#include "vtkImageChangeInformation.h"
#include "vtkOrientationMarkerWidget.h"
#include "vtkAnnotatedCubeActor.h"

```

```

#include "vtkAxesActor.h"
#include "vtkCaptionActor2D.h"
#include "vtkTextProperty.h"
#include "vtkPropAssembly.h"

#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkStringArray.h"

#include "gdcmSystem.h"
#include "gdcmDirectory.h"
#include "gdcmIPPSorter.h"

#ifdef vtkFloatingPointType
#define vtkFloatingPointType float
#endif

//-----
class vtkOrthoPlanesCallback : public vtkCommand
{
public:
    static vtkOrthoPlanesCallback *New()
    { return new vtkOrthoPlanesCallback; }

    void Execute( vtkObject *caller, unsigned long vtkNotUsed( event ),
                  void *callData )
    {
        vtkImagePlaneWidget* self =
            reinterpret_cast< vtkImagePlaneWidget* >( caller );
        if(!self) return;

        double* wl = static_cast<double*>( callData );

        if ( self == this->WidgetX )
        {
            this->WidgetY->SetWindowLevel(wl[0],wl[1],1);
            this->WidgetZ->SetWindowLevel(wl[0],wl[1],1);
        }
        else if( self == this->WidgetY )
        {
            this->WidgetX->SetWindowLevel(wl[0],wl[1],1);
            this->WidgetZ->SetWindowLevel(wl[0],wl[1],1);
        }
        else if (self == this->WidgetZ)
        {
            this->WidgetX->SetWindowLevel(wl[0],wl[1],1);
            this->WidgetY->SetWindowLevel(wl[0],wl[1],1);
        }
    }

    vtkOrthoPlanesCallback():WidgetX( 0 ), WidgetY( 0 ), WidgetZ ( 0 ) {}

    vtkImagePlaneWidget* WidgetX;
    vtkImagePlaneWidget* WidgetY;
    vtkImagePlaneWidget* WidgetZ;
};

int main( int argc, char *argv[] )
{
    //char* fname = vtkTestUtilities::ExpandDataFileName(argc, argv, "Data/headsq/quarter");

    //vtkVolume16Reader* v16 = vtkVolume16Reader::New();
    // v16->SetDataDimensions( 64, 64);
    // v16->SetDataByteOrderToLittleEndian();
    // v16->SetImageRange( 1, 93);
    // v16->SetDataSpacing( 3.2, 3.2, 1.5);
    // v16->SetFilePrefix( fname );
    // v16->SetDataMask( 0x7fff);
    // v16->Update();
    std::vector<std::string> filenames;
    if( argc < 2 )
    {
        std::cerr << argv[0] << " filename1.dcm [filename2.dcm ...]\n";
        return 1;
    }
    else
    {
        // Is it a single directory ? If so loop over all files contained in it:
        const char *filename = argv[1];
        if( argc == 2 && gdcm::System::FileIsDirectory( filename ) )
        {

```

```

std::cout << "Loading directory: " << filename << std::endl;
bool recursive = false;
gdcmm::Directory d;
d.Load(filename, recursive);
gdcmm::Directory::FileNamesType const &files = d.
GetFileNames();
for( gdcmm::Directory::FileNamesType::const_iterator it = files.begin(); it != files.end(); ++it )
{
    filenames.push_back( it->c_str() );
}
}
else // list of files passed directly on the cmd line:
    // discard non-existing or directory
    {
        for(int i=1; i < argc; ++i)
        {
            filename = argv[i];
            if( gdcmm::System::FileExists( filename ) )
            {
                if( gdcmm::System::FileIsDirectory( filename ) )
                {
                    std::cerr << "Discarding directory: " << filename << std::endl;
                }
                else
                {
                    filenames.push_back( filename );
                }
            }
            else
            {
                std::cerr << "Discarding non existing file: " << filename << std::endl;
            }
        }
    }
//names->Print( std::cout );
}

vtkGDCMImageReader * reader = vtkGDCMImageReader::New();
double ippzspacing;
if( filenames.size() > 1 )
{
    //gdcmm::Trace::DebugOn();
    //gdcmm::Trace::WarningOn();
    gdcmm::IPPSorter s;
    s.SetComputeZSpacing( true );
    s.SetZSpacingTolerance( 1e-3 );
    bool b = s.Sort( filenames );
    if( !b )
    {
        std::cerr << "Failed to sort files" << std::endl;
        return 1;
    }
    std::cout << "Sorting succeeded:" << std::endl;
    s.Print( std::cout );

    std::cout << "Found z-spacing:" << std::endl;
    std::cout << s.GetZSpacing() << std::endl;
    ippzspacing = s.GetZSpacing();

    const std::vector<std::string> & sorted = s.GetFileNames();
    vtkStringArray *files = vtkStringArray::New();
    std::vector< std::string >::const_iterator it = sorted.begin();
    for( ; it != sorted.end(); ++it )
    {
        const std::string &f = *it;
        files->InsertNextValue( f.c_str() );
    }
    reader->SetFileNames( files );
    //reader->SetFileLowerLeft( 1 );
    reader->Update(); // important
    files->Delete();
}
else
{
    reader->SetFileName( argv[1] );
    reader->Update(); // important
    ippzspacing = reader->GetOutput()->GetSpacing()[2];
    ippzspacing = 4;
}

//reader->GetOutput()->Print( std::cout );

```

```

//vtkFloatingPointType range[2];
//reader->GetOutput()->GetScalarRange(range);
//std::cout << "Range: " << range[0] << " " << range[1] << std::endl;

const vtkFloatingPointType *spacing = reader->GetOutput()->GetSpacing();

vtkImageChangeInformation *v16 = vtkImageChangeInformation::New();
#if (VTK_MAJOR_VERSION >= 6)
v16->SetInputConnection( reader->GetOutputPort() );
#else
v16->SetInput( reader->GetOutput() );
#endif
v16->SetOutputSpacing( spacing[0], spacing[1], ippszspacing );
v16->Update();

#if 0
vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
writer->SetInput( v16->GetOutput() );
writer->SetFileLowerLeft( reader->GetFileLowerLeft() );
writer->SetDirectionCosines( reader->GetDirectionCosines() );
writer->SetImageFormat( reader->GetImageFormat() );
writer->SetFileDimensionality( 3); //reader->GetFileDimensionality() );
writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
writer->SetShift( reader->GetShift() );
writer->SetScale( reader->GetScale() );
writer->SetFileName( "out.dcm" );
writer->Write();
#endif

vtkOutlineFilter* outline = vtkOutlineFilter::New();
outline->SetInputConnection(v16->GetOutputPort());

vtkPolyDataMapper* outlineMapper = vtkPolyDataMapper::New();
outlineMapper->SetInputConnection(outline->GetOutputPort());

vtkActor* outlineActor = vtkActor::New();
outlineActor->SetMapper( outlineMapper);

vtkRenderer* ren1 = vtkRenderer::New();
vtkRenderer* ren2 = vtkRenderer::New();

vtkRenderWindow* renWin = vtkRenderWindow::New();
renWin->AddRenderer(ren2);
renWin->AddRenderer(ren1);

vtkRenderWindowInteractor* iren = vtkRenderWindowInteractor::New();
iren->SetRenderWindow(renWin);

vtkCellPicker* picker = vtkCellPicker::New();
picker->SetTolerance(0.005);

vtkProperty* ipwProp = vtkProperty::New();
//assign default props to the ipw's texture plane actor

vtkImagePlaneWidget* planeWidgetX = vtkImagePlaneWidget::New();
planeWidgetX->SetInteractor( iren);
planeWidgetX->SetKeyPressActivationValue('x');
planeWidgetX->SetPicker(picker);
planeWidgetX->RestrictPlaneToVolumeOn();
planeWidgetX->GetPlaneProperty()->SetColor(1,0,0);
planeWidgetX->SetTexturePlaneProperty(ipwProp);
planeWidgetX->TextureInterpolateOff();
planeWidgetX->SetResliceInterpolateToNearestNeighbour();
#if (VTK_MAJOR_VERSION >= 6)
planeWidgetX->SetInputConnection(v16->GetOutputPort());
#else
planeWidgetX->SetInput(v16->GetOutput());
#endif
#endif
planeWidgetX->SetPlaneOrientationToXAxes();
//planeWidgetX->SetSliceIndex(32);
planeWidgetX->DisplayTextOn();
planeWidgetX->On();
planeWidgetX->InteractionOff();
planeWidgetX->InteractionOn();

vtkImagePlaneWidget* planeWidgetY = vtkImagePlaneWidget::New();
planeWidgetY->SetInteractor( iren);
planeWidgetY->SetKeyPressActivationValue('y');
planeWidgetY->SetPicker(picker);
planeWidgetY->GetPlaneProperty()->SetColor(1,1,0);

```

```

    planeWidgetY->SetTexturePlaneProperty(ipwProp);
    planeWidgetY->TextureInterpolateOn();
    planeWidgetY->SetResliceInterpolateToLinear();
    #if (VTK_MAJOR_VERSION >= 6)
        planeWidgetY->SetInputConnection(vl6->GetOutputPort());
    #else
        planeWidgetY->SetInput(vl6->GetOutput());
    #endif
    planeWidgetY->SetPlaneOrientationToYAxes();
    //planeWidgetY->SetSlicePosition(102.4);
    planeWidgetY->SetLookupTable(planeWidgetX->GetLookupTable());
    planeWidgetY->DisplayTextOn();
    planeWidgetY->UpdatePlacement();
    planeWidgetY->On();

    vtkImagePlaneWidget* planeWidgetZ = vtkImagePlaneWidget::New();
    planeWidgetZ->SetInteractor(iren);
    planeWidgetZ->SetKeyPressActivationValue('z');
    planeWidgetZ->SetPicker(picker);
    planeWidgetZ->GetPlaneProperty()->SetColor(0,0,1);
    planeWidgetZ->SetTexturePlaneProperty(ipwProp);
    planeWidgetZ->TextureInterpolateOn();
    planeWidgetZ->SetResliceInterpolateToCubic();
    #if (VTK_MAJOR_VERSION >= 6)
        planeWidgetZ->SetInputConnection(vl6->GetOutputPort());
    #else
        planeWidgetZ->SetInput(vl6->GetOutput());
    #endif
    planeWidgetZ->SetPlaneOrientationToZAxes();
    //planeWidgetZ->SetSliceIndex(25);
    planeWidgetZ->SetLookupTable(planeWidgetX->GetLookupTable());
    planeWidgetZ->DisplayTextOn();
    planeWidgetZ->On();

    vtkImageOrthoPlanes* orthoPlanes = vtkImageOrthoPlanes::New();
    orthoPlanes->SetPlane(0, planeWidgetX);
    orthoPlanes->SetPlane(1, planeWidgetY);
    orthoPlanes->SetPlane(2, planeWidgetZ);
    orthoPlanes->ResetPlanes();

    vtkOrthoPlanesCallback* cbk = vtkOrthoPlanesCallback::New();
    cbk->WidgetX = planeWidgetX;
    cbk->WidgetY = planeWidgetY;
    cbk->WidgetZ = planeWidgetZ;
    planeWidgetX->AddObserver(vtkCommand::EndWindowLevelEvent, cbk);
    planeWidgetY->AddObserver(vtkCommand::EndWindowLevelEvent, cbk);
    planeWidgetZ->AddObserver(vtkCommand::EndWindowLevelEvent, cbk);
    cbk->Delete();

    double wl[2];
    planeWidgetZ->GetWindowLevel(wl);

    // Add a 2D image to test the GetReslice method
    //
    vtkImageMapToColors* colorMap = vtkImageMapToColors::New();
    colorMap->PassAlphaToOutputOff();
    colorMap->SetActiveComponent(0);
    colorMap->SetOutputFormatToLuminance();
    #if (VTK_MAJOR_VERSION >= 6)
        colorMap->SetInputData(planeWidgetZ->GetResliceOutput());
    #else
        colorMap->SetInput(planeWidgetZ->GetResliceOutput());
    #endif
    colorMap->SetLookupTable(planeWidgetX->GetLookupTable());

    vtkImageActor* imageActor = vtkImageActor::New();
    imageActor->PickableOff();
    #if (VTK_MAJOR_VERSION >= 6)
        imageActor->SetInputData(colorMap->GetOutput());
    #else
        imageActor->SetInput(colorMap->GetOutput());
    #endif

    // Add the actors
    //
    ren1->AddActor(outlineActor);
    ren2->AddActor(imageActor);

    ren1->SetBackground(0.1, 0.1, 0.2);
    ren2->SetBackground(0.2, 0.1, 0.2);

```



```

renWin->SetSize( 600, 350);

ren1->SetViewport(0,0,0.58333,1);
ren2->SetViewport(0.58333,0,1,1);

// Set the actors' postions
//
renWin->Render();
//iren->SetEventPosition( 175,175);
//iren->SetKeyCode('r');
//iren->InvokeEvent(vtkCommand::CharEvent,NULL);
//iren->SetEventPosition( 475,175);
//iren->SetKeyCode('r');
//iren->InvokeEvent(vtkCommand::CharEvent,NULL);
//renWin->Render();

//ren1->GetActiveCamera()->Elevation(110);
//ren1->GetActiveCamera()->SetViewUp(0, 0, -1);
//ren1->GetActiveCamera()->Azimuth(45);
//ren1->GetActiveCamera()->Dolly(1.15);
ren1->ResetCameraClippingRange();

vtkAnnotatedCubeActor* cube = vtkAnnotatedCubeActor::New();
cube->SetXPlusFaceText ( "R" );
cube->SetXMinusFaceText ( "L" );
cube->SetYPlusFaceText ( "A" );
cube->SetYMinusFaceText ( "P" );
cube->SetZPlusFaceText ( "H" );
cube->SetZMinusFaceText ( "F" );
cube->SetFaceTextScale( 0.666667 );

vtkAxesActor* axes2 = vtkAxesActor::New();

vtkMatrix4x4 *invert = vtkMatrix4x4::New();
invert->DeepCopy( reader->GetDirectionCosines() );
invert->Invert();

// simulate a left-handed coordinate system
//
vtkTransform *transform = vtkTransform::New();
transform->Identity();
//transform->RotateY(90);
transform->Concatenate(invert);
axes2->SetShaftTypeToCylinder();
axes2->SetUserTransform( transform );
cube->GetAssembly()->SetUserTransform( transform );

axes2->SetTotalLength( 1.5, 1.5, 1.5 );
axes2->SetCylinderRadius( 0.500 * axes2->GetCylinderRadius() );
axes2->SetConeRadius ( 1.025 * axes2->GetConeRadius() );
axes2->SetSphereRadius ( 1.500 * axes2->GetSphereRadius() );

vtkTextProperty* tprop = axes2->GetXAxisCaptionActor2D()->
    GetCaptionTextProperty();
tprop->ItalicOn();
tprop->ShadowOn();
tprop->SetFontFamilyToTimes();

axes2->GetYAxisCaptionActor2D()->GetCaptionTextProperty()->ShallowCopy( tprop );
axes2->GetZAxisCaptionActor2D()->GetCaptionTextProperty()->ShallowCopy( tprop );

vtkPropAssembly* assembly = vtkPropAssembly::New();
assembly->AddPart( axes2 );
assembly->AddPart( cube );

vtkOrientationMarkerWidget* widget = vtkOrientationMarkerWidget::New();
widget->SetOutlineColor( 0.9300, 0.5700, 0.1300 );
widget->SetOrientationMarker( assembly );
widget->SetInteractor( iren );
widget->SetViewport( 0.0, 0.0, 0.4, 0.4 );
widget->SetEnabled( 1 );
widget->InteractiveOff();
widget->InteractiveOn();

// Playback recorded events
//
//vtkInteractorEventRecorder *recorder = vtkInteractorEventRecorder::New();
//recorder->SetInteractor(iren);
//recorder->ReadFromInputStringOn();
//recorder->SetInputString(IOEventLog);

```

```

// Interact with data
// Render the image
//
iren->Initialize();
renWin->Render();

// Test SetKeyPressActivationValue for one of the widgets
//
//iren->SetKeyCode('z');
//iren->InvokeEvent(vtkCommand::CharEvent,NULL);
//iren->SetKeyCode('z');
//iren->InvokeEvent(vtkCommand::CharEvent,NULL);

//int retVal = vtkRegressionTestImage( renWin );
//
//if ( retVal == vtkRegressionTester::DO_INTERACTOR)
//{
//    iren->Start();
//}

// Clean up
//
//recorder->Off();
//recorder->Delete();

ipwProp->Delete();
orthoPlanes->Delete();
planeWidgetX->Delete();
planeWidgetY->Delete();
planeWidgetZ->Delete();
colorMap->Delete();
imageActor->Delete();
picker->Delete();
outlineActor->Delete();
outlineMapper->Delete();
outline->Delete();
iren->Delete();
renWin->Delete();
ren1->Delete();
ren2->Delete();
v16->Delete();
reader->Delete();

return 0;
}

```

## 29.65 gdcmmreslice.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"

#include "vtkRenderer.h"
#include "vtkAssembly.h"
#include "vtkImageFlip.h"
#include "vtkImageReslice.h"
#include "vtkRenderWindow.h"
#include "vtkAnnotatedCubeActor.h"
#include "vtkTransform.h"
#include "vtkAxesActor.h"
#include "vtkTextProperty.h"
#include "vtkCaptionActor2D.h"
#include "vtkPropAssembly.h"
#include "vtkOrientationMarkerWidget.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkPolyDataMapper.h"

```

```

#include "vtkActor.h"
#include "vtkImageData.h"
#include "vtkLookupTable.h"
#include "vtkTexture.h"
#include "vtkPlaneSource.h"

int main( int argc, char *argv[] )
{
    if( argc < 2 ) return 1;
    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( argv[1] );
    //reader->FileLowerLeftOn();
    reader->Update();

    vtkImageFlip *flip = vtkImageFlip::New();
    #if (VTK_MAJOR_VERSION >= 6)
        flip->SetInputConnection(reader->GetOutputPort());
    #else
        flip->SetInput(reader->GetOutput());
    #endif
    flip->SetFilteredAxis(0);
    flip->Update();

    vtkImageReslice *reslice = vtkImageReslice::New();
    //reslice->SetInput(reader->GetOutput());
    #if (VTK_MAJOR_VERSION >= 6)
        reslice->SetInputConnection(flip->GetOutputPort());
    #else
        reslice->SetInput(flip->GetOutput());
    #endif
    //reslice->SetResliceAxesDirectionCosines()
    reader->GetDirectionCosines()->Print(std::cout);
    vtkMatrix4x4 *invert = vtkMatrix4x4::New();
    invert->DeepCopy( reader->GetDirectionCosines() );
    invert->Invert();

    //reslice->SetResliceAxes( reader->GetDirectionCosines() );
    reslice->SetResliceAxes( invert );
    reslice->Update();
    vtkImageData* ima = reslice->GetOutput();

    vtkLookupTable* table = vtkLookupTable::New();
    table->SetNumberOfColors(1000);
    table->SetTableRange(0,1000);
    table->SetSaturationRange(0,0);
    table->SetHueRange(0,1);
    table->SetValueRange(0,1);
    table->SetAlphaRange(1,1);
    table->Build();

    // Texture
    vtkTexture* texture = vtkTexture::New();
    #if (VTK_MAJOR_VERSION >= 6)
        texture->SetInputData(ima);
    #else
        texture->SetInput(ima);
    #endif
    texture->InterpolateOn();
    texture->SetLookupTable(table);

    // PlaneSource
    vtkPlaneSource* plane = vtkPlaneSource::New();

    // PolyDataMapper
    vtkPolyDataMapper *planeMapper = vtkPolyDataMapper::New();
    #if (VTK_MAJOR_VERSION >= 6)
        planeMapper->SetInputConnection(plane->GetOutputPort());
    #else
        planeMapper->SetInput(plane->GetOutput());
    #endif

    // Actor
    vtkActor* planeActor = vtkActor::New();
    planeActor->SetTexture(texture);
    planeActor->SetMapper(planeMapper);
    planeActor->PickableOn();

    // Final rendering with simple interactor:
    vtkRenderer *ren = vtkRenderer::New();
    vtkRenderWindow *renwin = vtkRenderWindow::New();
    renwin->AddRenderer(ren);

```

```

vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
iren->SetRenderWindow(renwin);
ren->AddActor(planeActor);
ren->SetBackground(0,0,0.5);

// DICOM is RAH:
vtkAnnotatedCubeActor* cube = vtkAnnotatedCubeActor::New();
cube->SetXPlusFaceText ( "R" );
cube->SetXMinusFaceText ( "L" );
cube->SetYPlusFaceText ( "A" );
cube->SetYMinusFaceText ( "P" );
cube->SetZPlusFaceText ( "H" );
cube->SetZMinusFaceText ( "F" );

vtkAxesActor* axes2 = vtkAxesActor::New();

vtkTransform *transform = vtkTransform::New();
transform->Identity();
//reader->GetDirectionCosines()->Print(std::cout);
transform->Concatenate(invert);
//axes2->SetShaftTypeToCylinder();
axes2->SetUserTransform( transform );
cube->GetAssembly()->SetUserTransform( transform ); // cant get it to work

vtkPropAssembly* assembly = vtkPropAssembly::New();
assembly->AddPart( axes2 );
assembly->AddPart( cube );

vtkOrientationMarkerWidget* widget = vtkOrientationMarkerWidget::New();
widget->SetOrientationMarker( assembly );
widget->SetInteractor( iren );
widget->SetEnabled( 1 );
widget->InteractiveOff();
widget->InteractiveOn();

renwin->Render();
iren->Start();

// Clean up:
reader->Delete();
table->Delete();
texture->Delete();
plane->Delete();
planeMapper->Delete();
planeActor->Delete();
ren->Delete();
renwin->Delete();
iren->Delete();

return 0;
}

```

## 29.66 gdcmrtnonplan.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcms.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkImageData.h"
#include "vtkPointData.h"
#include "vtkPolyData.h"
#include "vtkProperty.h"
#include "vtkPolyDataMapper.h"
#include "vtkActor.h"
#include "vtkRenderer.h"
#include "vtkCellArray.h"
#include "vtkPoints.h"
#include "vtkDoubleArray.h"

```

```

#include <vtkXMLImageDataWriter.h>
#include <vtkXMLPolyDataWriter.h>
#include <vtkRenderWindowInteractor.h>
#include <vtkImageColorViewer.h>

#include "gdcmmReader.h"
#include "gdcmmAttribute.h"

/*
This example is just for fun. We found a RT Ion Plan Storage and simply extracted the viz stuff for VTK

RTIonPlanStorage, // 1.2.840.10008.5.1.4.1.1.481.8
*/
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " filename.dcm outfile.vti\n";
        return 1;
    }
    const char * filename = argv[1];
    const char * outfilename = argv[2];
    const char * outfilename2 = argv[3];

    gdcmm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcmm::MediaStorage ms;
    ms.SetFromFile( reader.GetFile() );
    if( ms != gdcmm::MediaStorage::RTIonPlanStorage )
    {
        return 1;
    }

    /*
(300a,03a2) SQ # u/l,1 Ion Beam Sequence
(ffff,e000) na (Item with undefined length)
(0008,1040) LO [Test] # 4,1 Institutional Department Name
(300a,00b2) SH (no value) # 0,1 Treatment Machine Name
(300a,00b3) CS [MU] # 2,1 Primary Dosimeter Unit
(300a,00c0) IS [1 ] # 2,1 Beam Number
(300a,00c2) LO [1 ] # 2,1 Beam Name
(300a,00c4) CS [STATIC] # 6,1 Beam Type
(300a,00c6) CS [PROTON] # 6,1 Radiation Type
(300a,00ce) CS [TREATMENT ] # 10,1 Treatment Delivery Type
(300a,00d0) IS [0 ] # 2,1 Number of Wedges
(300a,00e0) IS [1 ] # 2,1 Number of Compensators
(300a,00ed) IS [0 ] # 2,1 Number of Boli
(300a,00f0) IS [1 ] # 2,1 Number of Blocks
(300a,0110) IS [2 ] # 2,1 Number of Control Points
(300a,02ea) SQ # u/l,1 Ion Range Compensator Sequence
(ffff,e000) na (Item with undefined length)
(300a,00e1) SH [lucite] # 6,1 Material ID
(300a,00e4) IS [1 ] # 2,1 Compensator Number
(300a,00e5) SH [75hdhe5 ] # 8,1 Compensator ID
(300a,00e7) IS [35] # 2,1 Compensator Rows
(300a,00e8) IS [37] # 2,1 Compensator Columns
(300a,00e9) DS [3.679991\4.249288 ] # 18,2 Compensator Pixel Spacing
(300a,00ea) DS [-76.00\62.50] # 12,2 Compensator Position
(300a,00ec) DS
[52.13\52.13\52.13\53.18\54.04\54.04\47.11\40.06\40.06\38.79\34.87\33.28\33.28\
33.28\35.43\35.43\34.54\34.54\34.71\36.10\38.62\44.88\44.88\44.88\45.00\45.00\45.00\45.66\45.66\46.42\39.77\39.77\39.77\39.77\
Data
(300a,02e0) CS [ABSENT] # 6,1 Compensator Divergence
(300a,02e1) CS [SOURCE_SIDE ] # 12,1 Compensator Mounting Position
(300a,02e4) FL 39.2 # 4,1 Isocenter to Compensator Tray
Distance
(300a,02e5) FL 2.12 # 4,1 Compensator Column Offset
(300a,02e8) FL 4.76 # 4,1 Compensator Milling Tool Diameter
(ffff,e00d)
*/
    const gdcmm::DataSet& ds = reader.GetFile().GetDataSet();
    gdcmm::Tag tbeamsq(0x300a,0x03a2);
    if( !ds.FindDataElement( tbeamsq ) )
    {
        return 1;
    }
}

```

```

const gdcm::DataElement &beamsq = ds.GetDataElement( tbeamsq );
//std::cout << beamsq << std::endl;
gdcm::SmartPointer<gdcm::SequenceOfItems> sqi = beamsq.
    GetValueAsSQ();
if( !sqi || !sqi->GetNumberOfItems() )
{
    return 1;
}

//for(unsigned int pd = 0; pd < sqi->GetNumberOfItems(); ++pd)
// {
    //const gdcm::Item & item = sqi->GetItem(1); // Item start at #1
    const gdcm::Item & item = sqi->GetItem(1); // Item start at #1
    const gdcm::DataSet& nestedds = item.GetNestedDataSet();
    //std::cout << nestedds << std::endl;
    gdcm::Tag tcompensatorsq(0x300a,0x02ea);
    if( !nestedds.FindDataElement( tcompensatorsq ) )
    {
        return 1;
    }
    const gdcm::DataElement &compensatorsq = nestedds.
        GetDataElement( tcompensatorsq );
    //std::cout << compensatorsq << std::endl;
    gdcm::SmartPointer<gdcm::SequenceOfItems> ssqi = compensatorsq.
        GetValueAsSQ();
    const gdcm::Item & item2 = ssqi->GetItem(1); // Item start at #1
    const gdcm::DataSet& nestedds2 = item2.GetNestedDataSet();
    //std::cout << nestedds2 << std::endl;
    gdcm::Tag tcompensatorthicknessdata(0x300a,0x00ec);
    if( !nestedds2.FindDataElement( tcompensatorthicknessdata ) )
    {
        return 1;
    }
    const gdcm::DataElement &compensatorthicknessdata = nestedds2.
        GetDataElement( tcompensatorthicknessdata );
    // std::cout << compensatorthicknessdata << std::endl;
    gdcm::Attribute<0x300a,0x00ec> at;
    at.SetFromDataElement( compensatorthicknessdata );
    const double* pts = at.GetValues();
    //      (300a,00e7) IS [35]                                # 2,1 Compensator Rows
    gdcm::Attribute<0x300a,0x00e7> at1;
    const gdcm::DataElement &compensatorrows = nestedds2.
        GetDataElement( at1.GetTag() );
    at1.SetFromDataElement( compensatorrows );
    std::cout << at1.GetValue() << std::endl;
    //      (300a,00e8) IS [37]                                # 2,1 Compensator Columns
    gdcm::Attribute<0x300a,0x00e8> at2;
    const gdcm::DataElement &compensatorcols = nestedds2.
        GetDataElement( at2.GetTag() );
    at2.SetFromDataElement( compensatorcols );
    std::cout << at2.GetValue() << std::endl;

    // (300a,00e9) DS [3.679991\4.249288 ]                    # 18,2 Compensator Pixel Spacing
    gdcm::Attribute<0x300a,0x00e9> at3;
    const gdcm::DataElement &compensatorpixelspacing = nestedds2.
        GetDataElement( at3.GetTag() );
    at3.SetFromDataElement( compensatorpixelspacing );
    std::cout << at3.GetValue(0) << std::endl;
    // (300a,00ea) DS [-76.00\62.50]                            # 12,2 Compensator Position
    gdcm::Attribute<0x300a,0x00ea> at4;
    const gdcm::DataElement &compensatorposition = nestedds2.
        GetDataElement( at4.GetTag() );
    at4.SetFromDataElement( compensatorposition );
    std::cout << at4.GetValue(0) << std::endl;

    vtkDoubleArray *d = vtkDoubleArray::New();
    d->SetArray( (double*)pts , at1.GetValue() * at2.GetValue() , 0 );

    vtkImageData *img = vtkImageData::New();
    img->Initialize();
    img->SetDimensions( at2.GetValue(), at1.GetValue(), 1 );
    //imgb->SetExtent(1, xdim, 1, ydim, 1, zdim);
    #if (VTK_MAJOR_VERSION >= 6)
        assert(0);
    #else
        img->SetScalarTypeToDouble();
    #endif
    img->SetSpacing( at3.GetValue(1), at3.GetValue(0), 1); // FIXME image is upside down
    img->SetOrigin( at4.GetValue(0), at4.GetValue(1), 1);
    #if (VTK_MAJOR_VERSION >= 6)
        assert(0);
    #endif

```

```

#else
    img->SetNumberOfScalarComponents(1);
#endif
img->GetPointData()->SetScalars(d);

#if (VTK_MAJOR_VERSION >= 6)
#else
    img->Update();
#endif
img->Print(std::cout);

    vtkXMLImageDataWriter *writeb= vtkXMLImageDataWriter::New();
#if (VTK_MAJOR_VERSION >= 6)
    writeb->SetInputData( img );
#else
    writeb->SetInput( img );
#endif
writeb->SetFileName( outfilename );
writeb->Write();

/*
    (300a,03a6) SQ                                     # u/1,1 Ion Block Sequence
    (fffe,e000) na (Item with undefined length)
    (300a,00e1) SH [brass ]                             # 6,1 Material ID
    (300a,00f7) FL 95.03                               # 4,1 Isocenter to Block Tray Distance
    (300a,00f8) CS [APERTURE]                           # 8,1 Block Type
    (300a,00fa) CS [ABSENT]                             # 6,1 Block Divergence
    (300a,00fb) CS [SOURCE_SIDE ]                       # 12,1 Block Mounting Position
    (300a,00fc) IS [1 ]                                 # 2,1 Block Number
    (300a,0100) DS [50.00 ]                             # 6,1 Block Thickness
    (300a,0104) IS [179 ]                               # 4,1 Block Number of Points
    (300a,0106) DS
    [1.7\50.0\14.3\50.0\16.7\49.4\18.7\48.2\19.4\47.7\20.1\47.1\21.0\47.0\22.3\47.0\23.7\
    46.8\25.7\46.2\27.0\45.6\27.2\45.4\28.2\44.6\28.9\44.2\29.7\43.9\31.5\43.5\33.0\42.8\33.7\42.4\35.2\41.3\38.2\40.4\39.6\39.7\
    2\37.4\43.0\37.1\44.7\36] # 1934,2-2n Block Data
    (fffe,e00d)
    (fffe,e0dd)

*/

gdcm::Tag tblocksq(0x300a,0x03a6);
if( !nestedds.FindDataElement( tblocksq ) )
{
    return 1;
}
const gdcm::DataElement &tblocksq = nestedds.GetDataElement( tblocksq );
//std::cout << tblocksq << std::endl;
gdcm::SmartPointer<gdcm::SequenceOfItems> sssqi = tblocksq.
    GetValueAsSQ();
const gdcm::Item &item3 = sssqi->GetItem(1); // Item start at #1
const gdcm::DataSet &nestedds3 = item3.GetNestedDataSet();

gdcm::Tag tblockdata(0x300a,0x0106);
if( !nestedds3.FindDataElement( tblockdata ) )
{
    return 1;
}
const gdcm::DataElement &tblockdata = nestedds3.
    GetDataElement( tblockdata );
// std::cout << tblockdata << std::endl;
gdcm::Attribute<0x300a,0x0106> at_;
at_.SetFromDataElement( tblockdata );

vtkDoubleArray *scalars = vtkDoubleArray::New();
scalars->SetNumberOfComponents(3);

gdcm::Attribute<0x300a,0x0104> bnpts; // IS [179 ]
    # 4,1 Block Number of Points
if( !nestedds3.FindDataElement( bnpts.GetTag() ) )
{
    return 1;
}
const gdcm::DataElement &tblocknpts = nestedds3.
    GetDataElement( bnpts.GetTag() );
bnpts.SetFromDataElement( tblocknpts );
//std::cout << bnpts.GetValue() << std::endl;

vtkPolyData *output = vtkPolyData::New();
vtkPoints *newPts = vtkPoints::New();
vtkCellArray *polys = vtkCellArray::New();
const double *ptr = at_.GetValues();
//unsigned int npts = bnpts.GetNumberOfValues() / 2;
unsigned int npts = bnpts.GetValue();

```

```

vtkIdType *ptIds = new vtkIdType[npts];
for(unsigned int i = 0; i < npts; ++i)
{
    float x[3] = {};
    x[0] = (float)ptr[2*i+0];
    x[1] = (float)ptr[2*i+1];
    //x[2] = pts[i+2];
    vtkIdType ptId = newPts->InsertNextPoint( x );
    //std::cout << x[0] << ", " << x[1] << ", " << x[2] << std::endl;
    ptIds[i] = ptId;
}
vtkIdType cellId = polys->InsertNextCell(npts , ptIds);
(void)cellId;
delete[] ptIds;

output->SetPoints(newPts);
newPts->Delete();
output->SetPolys(polys);
polys->Delete();
//output->GetCellData()->SetScalars(scalars);
//scalars->Delete();
#if (VTK_MAJOR_VERSION >= 6)
#else
    output->Update();
#endif
output->Print( std::cout );

// }

vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();

vtkImageColorViewer *viewer = vtkImageColorViewer::New();
#if (VTK_MAJOR_VERSION >= 6)
viewer->SetInputData(img);
#else
viewer->SetInput(img);
#endif
viewer->SetupInteractor(iren);
viewer->SetSize(600, 600);
viewer->GetRenderer()->ResetCameraClippingRange();
viewer->Render();
viewer->GetRenderer()->ResetCameraClippingRange();

vtkPolyDataMapper *cubeMapper = vtkPolyDataMapper::New();
//vtkPolyDataMapper2D* cubeMapper = vtkPolyDataMapper2D::New();
#if (VTK_MAJOR_VERSION >= 6)
cubeMapper->SetInputData( output );
#else
cubeMapper->SetInput( output );
#endif
cubeMapper->SetScalarRange(0,7);
vtkActor *cubeActor = vtkActor::New();
//vtkActor2D* cubeActor = vtkActor2D::New();
cubeActor->SetMapper(cubeMapper);
vtkProperty * property = cubeActor->GetProperty();
property->SetRepresentationToWireframe();

viewer->GetRenderer()->AddActor( cubeActor );

vtkXMLPolyDataWriter *writec= vtkXMLPolyDataWriter::New();
#if (VTK_MAJOR_VERSION >= 6)
writec->SetInputData( output );
#else
writec->SetInput( output );
#endif
writec->SetFileName( outfilename2 );
writec->Write();

iren->Initialize();
iren->Start();

return 0;
}

```



## 29.67 gdcmrtplan.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkImageData.h"
#include "vtkPointData.h"
#include "vtkPolyData.h"
#include "vtkProperty.h"
#include "vtkPolyDataMapper.h"
#include "vtkActor.h"
#include "vtkRenderer.h"
#include "vtkCellArray.h"
#include "vtkPoints.h"
#include "vtkDoubleArray.h"
#include <vtkXMLImageDataWriter.h>
#include <vtkRenderWindowInteractor.h>
#include <vtkImageColorViewer.h>

#include "gdcmReader.h"
#include "gdcmAttribute.h"

/*
This example is just for fun. We found a fake RT Ion Plan Storage and simply extracted the viz stuff for
  VTK
but this is rather a RT Plan storage
*/
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " filename.dcm outfile.vti\n";
        return 1;
    }
    const char * filename = argv[1];
    const char * outfilename = argv[2];

    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::MediaStorage ms;
    ms.SetFromFile( reader.GetFile() );
    if( ms != gdcm::MediaStorage::RTIonPlanStorage )
    {
        return 1;
    }

    /*
(300a,00b0) SQ                                     # u/1,1 Beam Sequence
  (fffe,e000) na (Item with undefined length)
    (300a,00b2) SH (no value)                       # 0,1 Treatment Machine Name
    (300a,00c0) IS [1 ]                             # 2,1 Beam Number
    (300a,00c2) LO [1 ]                             # 2,1 Beam Name
    (300a,00c4) CS [STATIC]                         # 6,1 Beam Type
    (300a,00c6) CS [PROTON]                         # 6,1 Radiation Type
    (300a,00ce) CS [TREATMENT ]                    # 10,1 Treatment Delivery Type
    (300a,00e0) IS [1 ]                             # 2,1 Number of Compensators
    (300a,00e3) SQ                                   # u/1,1 Compensator Sequence
      (fffe,e000) na (Item with undefined length)
        (300a,00e1) SH [lucite]                     # 6,1 Material ID
        (300a,00e4) IS [1 ]                         # 2,1 Compensator Number
        (300a,00e5) SH [75hdhe5 ]                   # 8,1 Compensator ID
        (300a,00e7) IS [35]                         # 2,1 Compensator Rows
        (300a,00e8) IS [37]                         # 2,1 Compensator Columns
        (300a,00e9) DS [3.679991\4.249288 ]         # 18,2 Compensator Pixel Spacing
        (300a,00ea) DS [-76.00\62.50]               # 12,2 Compensator Position
    */
}

```



```

vtkDoubleArray *d = vtkDoubleArray::New();
d->SetArray( (double*)pts , at1.GetValue() * at2.GetValue() , 0 );

vtkImageData *img = vtkImageData::New();
img->Initialize();
img->SetDimensions( at2.GetValue(), at1.GetValue(), 1 );
//img->SetExtent(1, xdim, 1, ydim, 1, zdim);
#if (VTK_MAJOR_VERSION >= 6)
    assert(0);
#else
    img->SetScalarTypeToDouble();
#endif
img->SetSpacing( at3.GetValue(1), at3.GetValue(0), 1); // FIXME image is upside down
img->SetOrigin( at4.GetValue(0), at4.GetValue(1), 1);
#if (VTK_MAJOR_VERSION >= 6)
    assert(0);
#else
    img->SetNumberOfScalarComponents(1);
#endif
img->GetPointData()->SetScalars(d);

vtkXMLImageDataWriter *writeb= vtkXMLImageDataWriter::New();
#if (VTK_MAJOR_VERSION >= 6)
    writeb->SetInputData( img );
#else
    writeb->SetInput( img );
#endif
writeb->SetFileName( outfilename );
writeb->Write();
/*
(300a,00f4) SQ # u/1,1 Block Sequence
    (fffe,e000) na (Item with undefined length)
    (300a,00e1) SH [brass ] # 6,1 Material ID
    (300a,00f8) CS [APERTURE] # 8,1 Block Type
    (300a,00fa) CS [ABSENT] # 6,1 Block Divergence
    (300a,00fb) CS [SOURCE_SIDE ] # 12,1 Block Mounting Position
    (300a,00fc) IS [1 ] # 2,1 Block Number
    (300a,0100) DS [50.00 ] # 6,1 Block Thickness
    (300a,0104) IS [179 ] # 4,1 Block Number of Points
    (300a,0106) DS
    [1.7\50.0\14.3\50.0\16.7\49.4\18.7\48.2\19.4\47.7\20.1\47.1\21.0\47.0\22.3\47.0\23.7\
    46.8\25.7\46.2\27.0\45.6\27.2\45.4\28.2\44.6\28.9\44.2\29.7\43.9\31.5\43.5\33.0\42.8\33.7\42.4\35.2\41.3\38.2\40.4\39.6\39.7\
    (fffe,e00d)
    (fffe,e000) na (Item with undefined length)
    (fffe,e00d)
    (fffe,e0dd)
*/
gdcmm::Tag tblocksq(0x300a,0x00f4);
if( !nestedds.FindDataElement( tblocksq ) )
{
    return 1;
}
const gdcmm::DataElement &blocksq = nestedds.GetDataElement( tblocksq );
//std::cout << blocksq << std::endl;
gdcmm::SmartPointer<gdcmm::SequenceOfItems> sssqi = blocksq.
    GetValueAsSQ();
const gdcmm::Item & item3 = sssqi->GetItem(1); // Item start at #1
const gdcmm::DataSet& nestedds3 = item3.GetNestedDataSet();

gdcmm::Tag tblockdata(0x300a,0x0106);
if( !nestedds3.FindDataElement( tblockdata ) )
{
    return 1;
}
const gdcmm::DataElement &tblockdata = nestedds3.
    GetDataElement( tblockdata );
// std::cout << tblockdata << std::endl;
gdcmm::Attribute<0x300a,0x0106> at_;
at_.SetFromDataElement( tblockdata );

vtkDoubleArray *scalars = vtkDoubleArray::New();
scalars->SetNumberOfComponents(3);

gdcmm::Attribute<0x300a,0x0104> bnpts; // IS [179 ] # 4,1 Block Number of
    Points
if( !nestedds3.FindDataElement( bnpts.GetTag() ) )
{
    return 1;
}
const gdcmm::DataElement &tblocknpts = nestedds3.

```

```

    GetDataElement( bnpts.GetTag() );
    bnpts.SetFromDataElement( blocknpts );
    std::cout << bnpts.GetValue() << std::endl;

    vtkPolyData *output = vtkPolyData::New();
    vtkPoints *newPts = vtkPoints::New();
    vtkCellArray *polys = vtkCellArray::New();
    const double *ptr = at_.GetValues();
    //unsigned int npts = bnpts.GetNumberOfValues() / 2;
    unsigned int npts = bnpts.GetValue();
    vtkIdType *ptIds = new vtkIdType[npts];
    for(unsigned int i = 0; i < npts; ++i)
    {
        float x[3] = {};
        x[0] = (float)ptr[2*i+0];
        x[1] = (float)ptr[2*i+1];
        //x[2] = ptr[i+2];
        vtkIdType ptId = newPts->InsertNextPoint( x );
        //std::cout << x[0] << ", " << x[1] << ", " << x[2] << std::endl;
        ptIds[i] = ptId;
    }
    vtkIdType cellId = polys->InsertNextCell(npts, ptIds);
    (void)cellId;
    delete[] ptIds;

    output->SetPoints(newPts);
    newPts->Delete();
    output->SetPolys(polys);
    polys->Delete();
    //output->GetCellData()->SetScalars(scalars);
    //scalars->Delete();
    #if (VTK_MAJOR_VERSION >= 6)
    #else
        output->Update();
    #endif
    output->Print( std::cout );

    // }

    vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();

    vtkImageColorViewer *viewer = vtkImageColorViewer::New();
    #if (VTK_MAJOR_VERSION >= 6)
        viewer->SetInputData(img);
    #else
        viewer->SetInput(img);
    #endif
    viewer->SetupInteractor(iren);
    viewer->SetSize(600, 600);
    viewer->Render();

    vtkPolyDataMapper *cubeMapper = vtkPolyDataMapper::New();
    //vtkPolyDataMapper2D* cubeMapper = vtkPolyDataMapper2D::New();
    #if (VTK_MAJOR_VERSION >= 6)
        cubeMapper->SetInputData( output );
    #else
        cubeMapper->SetInput( output );
    #endif
    cubeMapper->SetScalarRange(0,7);
    vtkActor *cubeActor = vtkActor::New();
    //vtkActor2D* cubeActor = vtkActor2D::New();
    cubeActor->SetMapper(cubeMapper);
    vtkProperty *property = cubeActor->GetProperty();
    property->SetRepresentationToWireframe();

    viewer->GetRenderer()->AddActor( cubeActor );

    iren->Initialize();
    iren->Start();

    return 0;
}

```

## 29.68 gdcmscene.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMPolyDataReader.h"
// #include "vtkGDCMPolyDataWriter.h"

#include "vtkAppendPolyData.h"
#include "vtkPolyDataWriter.h"
#include "vtkPolyDataMapper.h"
#include "vtkPolyDataMapper2D.h"
#include "vtkActor2D.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkRenderWindow.h"
#include "vtkRenderer.h"
#include "vtkCamera.h"
#include "vtkProperty.h"
#include "vtkProperty2D.h"

// gdcmDataExtra/gdcmNonImageData/exRT_Structure_Set_Storage.dcm
// gdcmDataExtra/gdcmNonImageData/RTSTRUCT_1.3.6.1.4.1.22213.1.1396.2.dcm
// gdcmDataExtra/gdcmNonImageData/RT/RTStruct.dcm

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " filename1.dcm\n";
        return 1;
    }
    const char * filename = argv[1];

    vtkGDCMPolyDataReader * reader =
        vtkGDCMPolyDataReader::New();
    reader->SetFileName( filename );
    reader->Update();

    // vtkGDCMPolyDataWriter * writer2 = vtkGDCMPolyDataWriter::New();
    // for(int num = 0; num < reader->GetNumberOfOutputPorts(); ++num )
    //     writer2->SetInput( num, reader->GetOutput(num) );
    // writer2->SetFileName( "rtstruct.dcm" );
    // writer2->Write();

    // print reader output:
    reader->Print( std::cout );
    // print first output:
    reader->GetOutput()->Print( std::cout );

    vtkAppendPolyData *append = vtkAppendPolyData::New();
    int n = reader->GetNumberOfOutputPorts();
    for(int i = 0; i < n; ++i)
    {
        #if (VTK_MAJOR_VERSION >= 6)
            append->AddInputConnection( reader->GetOutputPort(i) );
        #else
            append->AddInput( reader->GetOutput(i) );
        #endif
    }

    vtkPolyDataWriter * writer = vtkPolyDataWriter::New();
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputConnection( reader->GetOutputPort() );
    #else
        writer->SetInput( reader->GetOutput() );
    #endif
    writer->SetFileName( "rtstruct.vtk" );
    //writer->Write();
}

```

```

    // Now we'll look at it.
    vtkPolyDataMapper *cubeMapper = vtkPolyDataMapper::New();
    //vtkPolyDataMapper2D* cubeMapper = vtkPolyDataMapper2D::New();
    //cubeMapper->SetInput( reader->GetOutput() );
    if (VTK_MAJOR_VERSION >= 6)
        cubeMapper->SetInputConnection( append->GetOutputPort() );
    #else
        cubeMapper->SetInput( append->GetOutput() );
    #endif
    cubeMapper->SetScalarRange(0,7);
    vtkActor *cubeActor = vtkActor::New();
    //vtkActor2D* cubeActor = vtkActor2D::New();
    cubeActor->SetMapper(cubeMapper);
    vtkProperty *property = cubeActor->GetProperty();
    property->SetRepresentationToWireframe();
    //cubeActor->GetProperty()->SetColor(1, 0, 0);

    // The usual rendering stuff.
    // vtkCamera *camera = vtkCamera::New();
    // camera->SetPosition(1,1,1);
    // camera->SetFocalPoint(0,0,0);

    vtkRenderer *renderer = vtkRenderer::New();
    vtkRenderWindow *renWin = vtkRenderWindow::New();
    renWin->AddRenderer(renderer);

    vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
    iren->SetRenderWindow(renWin);

    renderer->AddActor(cubeActor);
    //renderer->AddActor2D(cubeActor);
    //renderer->SetActiveCamera(camera);
    renderer->ResetCamera();
    renderer->SetBackground(1,1,1);

    renWin->SetSize(300,300);

    // interact with data
    renWin->Render();
    iren->Start();

    reader->Delete();
    append->Delete();
    cubeMapper->Delete();
    cubeActor->Delete();
    // camera->Delete();
    renderer->Delete();
    renWin->Delete();
    iren->Delete();

    writer->Delete();

    return 0;
}

```

## 29.69 gdcmttexture.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"

#include "vtkRenderer.h"
#include "vtkAssembly.h"

```

```

#include "vtkRenderWindow.h"
#include "vtkAnnotatedCubeActor.h"
#include "vtkTransform.h"
#include "vtkAxesActor.h"
#include "vtkTextProperty.h"
#include "vtkCaptionActor2D.h"
#include "vtkPropAssembly.h"
#include "vtkOrientationMarkerWidget.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkPolyDataMapper.h"
#include "vtkActor.h"
#include "vtkImageData.h"
#include "vtkLookupTable.h"
#include "vtkTexture.h"
#include "vtkPlaneSource.h"

int main( int argc, char *argv[] )
{
    if( argc < 2 ) return 1;
    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( argv[1] );

    reader->Update();
    vtkImageData* ima = reader->GetOutput();

    vtkLookupTable* table = vtkLookupTable::New();
    table->SetNumberOfColors(1000);
    table->SetTableRange(0,1000);
    table->SetSaturationRange(0,0);
    table->SetHueRange(0,1);
    table->SetValueRange(0,1);
    table->SetAlphaRange(1,1);
    table->Build();

    // Texture
    vtkTexture* texture = vtkTexture::New();
    #if (VTK_MAJOR_VERSION >= 6)
        texture->SetInputData(ima);
    #else
        texture->SetInput(ima);
    #endif
    texture->InterpolateOn();
    texture->SetLookupTable(table);

    // PlaneSource
    vtkPlaneSource* plane = vtkPlaneSource::New();
    plane->SetOrigin( -0.5, -0.5, 0.0);
    plane->SetPoint1( 0.5, -0.5, 0.0);
    plane->SetPoint2( -0.5, 0.5, 0.0);

    // PolyDataMapper
    vtkPolyDataMapper *planeMapper = vtkPolyDataMapper::New();
    #if (VTK_MAJOR_VERSION >= 6)
        planeMapper->SetInputConnection(plane->GetOutputPort());
    #else
        planeMapper->SetInput(plane->GetOutput());
    #endif

    // Actor
    vtkActor* planeActor = vtkActor::New();
    planeActor->SetTexture(texture);
    planeActor->SetMapper(planeMapper);
    planeActor->PickableOn();

    // Final rendering with simple interactor:
    vtkRenderer *ren = vtkRenderer::New();
    vtkRenderWindow *renwin = vtkRenderWindow::New();
    renwin->AddRenderer(ren);
    vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
    iren->SetRenderWindow(renwin);
    ren->AddActor(planeActor);
    ren->SetBackground(0,0,0.5);

    vtkAnnotatedCubeActor* cube = vtkAnnotatedCubeActor::New();
    cube->SetXPlusFaceText ( "L" );
    cube->SetXMinusFaceText ( "R" );
    cube->SetYPlusFaceText ( "A" );
    cube->SetYMinusFaceText ( "P" );
    cube->SetZPlusFaceText ( "H" );
    cube->SetZMinusFaceText ( "F" );

```

```

vtkAxesActor* axes2 = vtkAxesActor::New();
// simulate a left-handed coordinate system
//
vtkTransform *transform = vtkTransform::New();
transform->Identity();
//transform->RotateY(180);
reader->GetDirectionCosines()->Print(std::cout);
transform->Concatenate(reader->GetDirectionCosines());
//axes2->SetShaftTypeToCylinder();
axes2->SetUserTransform( transform );
//cube->SetUserTransform( transform ); // cant get it to work
cube->GetAssembly()->SetUserTransform( transform ); // cant get it to work

vtkPropAssembly* assembly = vtkPropAssembly::New();
assembly->AddPart( axes2 );
assembly->AddPart( cube );

vtkOrientationMarkerWidget* widget = vtkOrientationMarkerWidget::New();
//widget->SetOutlineColor( 0.9300, 0.5700, 0.1300 );
widget->SetOrientationMarker( assembly );
widget->SetInteractor( iren );
//widget->SetViewport( 0.0, 0.0, 0.4, 0.4 );
widget->SetEnabled( 1 );
widget->InteractiveOff();
widget->InteractiveOn();

renwin->Render();
iren->Start();

// Clean up:
reader->Delete();
table->Delete();
texture->Delete();
plane->Delete();
planeMapper->Delete();
planeActor->Delete();
ren->Delete();
renwin->Delete();
iren->Delete();

return 0;
}

```

## 29.70 gdcmvolume.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkPiecewiseFunction.h"
#include "vtkColorTransferFunction.h"
#include "vtkVolume.h"
#include "vtkVolumeProperty.h"
#include "vtkVolumeTextureMapper3D.h"
#include "vtkFixedPointVolumeRayCastMapper.h"
#include "vtkInteractorStyleTrackballCamera.h"
#include "vtkRenderer.h"
#include "vtkRenderWindow.h"
#include "vtkImageClip.h"
#include "vtkRenderWindowInteractor.h"

// gdcmvolume gdcmData/GE_DLX-8-MONO2-Multiframe-Jpeg_Lossless.dcm
int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();

```



```

reader->SetFileName( argv[1] );
reader->Update();

// Create the renderers, render window, and interactor
vtkRenderWindow *renWin = vtkRenderWindow::New();
vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
iren->SetRenderWindow(renWin);
vtkRenderer *ren = vtkRenderer::New();
renWin->AddRenderer(ren);

// Create a transfer function mapping scalar value to opacity
vtkPiecewiseFunction *oTFun = vtkPiecewiseFunction::New();
//oTFun->AddSegment(0, 1.0, 256, 0.1);
oTFun->AddSegment(0, 1.0, 240, 0.1);

vtkColorTransferFunction *cTFun = vtkColorTransferFunction::New();
cTFun->AddRGBPoint( 0, 1.0, 1.0, 1.0 );
//cTFun->AddRGBPoint( 255, 1.0, 1.0, 1.0 );
cTFun->AddRGBPoint( 240, 1.0, 1.0, 1.0 );

// Need to crop to actually see minimum intensity
vtkImageClip *clip = vtkImageClip::New();
clip->SetInputConnection( reader->GetOutputPort() );
clip->SetOutputWholeExtent(0,66,0,66,30,37);
clip->ClipDataOn();

vtkVolumeProperty *property = vtkVolumeProperty::New();
property->SetScalarOpacity(oTFun);
property->SetColor(cTFun);
property->SetInterpolationTypeToLinear();

vtkFixedPointVolumeRayCastMapper *mapper = vtkFixedPointVolumeRayCastMapper::New();
mapper->SetBlendModeToMinimumIntensity();
mapper->SetInputConnection( reader->GetOutputPort() );

vtkVolume *volume = vtkVolume::New();
volume->SetMapper(mapper);
volume->SetProperty(property);

ren->AddViewProp(volume);

renWin->Render();
{
    iren->Start();
}

volume->Delete();
mapper->Delete();
property->Delete();
clip->Delete();
cTFun->Delete();
oTFun->Delete();
reader->Delete();
renWin->Delete();
iren->Delete();
ren->Delete();

return 0;
}

```

## 29.71 GenAIIVR.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

```

```

#include "gdcmReader.h"
#include "gdcmGlobal.h"
#include "gdcmDummyValueGenerator.h"
#include "gdcmMediaStorage.h"
#include "gdcmWriter.h"
#include "gdcmItem.h"
#include "gdcmImageReader.h"
#include "gdcmSequenceOfItems.h"
#include "gdcmFile.h"
#include "gdcmTag.h"
#include "gdcmDict.h"
#include "gdcmDictEntry.h"
#include "gdcmDicts.h"
#include "gdcmTransferSyntax.h"
#include "gdcmUIDGenerator.h"
#include "gdcmFileExplicitFilter.h"

#include <cstdlib>
#include <cstring>

gdcm::Tag FindTagFromVR(gdcm::Dict const &dict, gdcm::VR const &vr)
{
    using gdcm::Dict;
    Dict::ConstIterator beg = dict.Begin();
    Dict::ConstIterator end = dict.End();
    Dict::ConstIterator it;
    for( it = beg; it != end; ++it)
    {
        const gdcm::Tag &t = it->first;
        const gdcm::DictEntry &de = it->second;
        const gdcm::VR &vr_de = de.GetVR();
        if( vr == vr_de && !de.GetRetired() && t.GetGroup() >= 0x8 )
        {
            return t;
        }
    }
    return gdcm::Tag(0xffff,0xffff);
}

struct rnd_gen {
    rnd_gen(char const* r = "abcdefghijklmnopqrstuvwxy0123456789")
        : range(r), len(std::strlen(r)) { }

    char operator ()() const {
        return range[static_cast<std::size_t>(std::rand() * (1.0 / ((double)RAND_MAX + 1.0 )) * (double)len)];
    }
private:
    char const* range;
    std::size_t len;
};

/*
*/
int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " output.dcm" << std::endl;
        return 1;
    }
    const char *outfilename = argv[1];
    static const gdcm::Global &g = gdcm::Global::GetInstance();
    static const gdcm::Dicts &dicts = g.GetDicts();
    static const gdcm::Dict &pubdict = dicts.GetPublicDict();
    using gdcm::VR;
    using gdcm::Tag;

    gdcm::Writer w;

    gdcm::File &f = w.GetFile();
    gdcm::DataSet &ds = f.GetDataSet();

    gdcm::FileExplicitFilter fef;
    //fef.SetChangePrivateTags( true );
    fef.SetFile( w.GetFile() );
    if( !fef.Change() )
    {
        std::cerr << "Failed to change" << std::endl;
        return 1;
    }
}

```

```

gdcmm::SmartPointer<gdcmm::SequenceOfItems> sq = new
    gdcmm::SequenceOfItems();
sq->SetLengthToUndefined();

// gdcmm::DummyValueGenerator dvg;

const std::size_t len = 10;
char ss[len+1];
ss[len] = '\\0';

const char owner_str[] = "GDCM CONFORMANCE TESTS";
gdcmm::DataElement owner( gdcmm::Tag(0x4d4d, 0x10) );
owner.SetByteValue(owner_str, (uint32_t)strlen(owner_str));
owner.SetVR( gdcmm::VR::LO );

// Create an item
gdcmm::Item it;
it.SetVLToUndefined();
gdcmm::DataSet &nds = it.GetNestedDataSet();
// nds.Insert(owner);
// nds.Insert(de);

// Insert sequence into data set
gdcmm::DataElement des( gdcmm::Tag(0x4d4d,0x1001) );
des.SetVR(gdcmm::VR::SQ);
des.SetValue(*sq);
des.SetVLToUndefined();

ds.Insert(owner);
ds.Insert(des);

// avoid INVALID = 0
for(int i = 1; i < 27; ++i)
{
    VR vr = (VR::VRType)(1 << i);
    Tag t = FindTagFromVR( pubdict, vr );
    if( vr != VR::UN && vr != VR::SQ )
    {
        assert( t != Tag(0xffff,0xffff) );
        gdcmm::DataElement de( t );
        std::generate_n(ss, len, rnd_gen());
        de.SetVR( vr );
        de.SetByteValue( ss, (uint32_t)std::strlen( ss ) );
        nds.Insert( de );
    }
}
sq->AddItem(it);

// Make sure to override any UID stuff
gdcmm::UIDGenerator uid;
gdcmm::DataElement de( Tag(0x8,0x18) ); // SOP Instance UID
de.SetVR( VR::UI );
const char *u = uid.Generate();
de.SetByteValue( u, (uint32_t)strlen(u) );
ds.Insert( de );

de.SetTag( Tag(0x8,0x16) ); // SOP Class UID
de.SetVR( VR::UI );
gdcmm::MediaStorage ms( gdcmm::MediaStorage::RawDataStorage
    );
de.SetByteValue( ms.GetString(), (uint32_t)strlen(ms.GetString()));
ds.Insert( de );

gdcmm::FileMetaInformation &fmi = f.GetHeader();
//fmi.SetDataSetTransferSyntax( gdcmm::TransferSyntax::ImplicitVRLittleEndian );
fmi.SetDataSetTransferSyntax(
    gdcmm::TransferSyntax::ExplicitVRLittleEndian );

w.SetCheckFileMetaInformation( true );
w.SetFileName( outfilename );
if (!w.Write() )
{
    return 1;
}

return 0;
}

```

## 29.72 GenerateDICOMDIR.cs

This is a C# example on how to use DICOMDIRGenerator

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Simple C# example to show how to use DICOMDIRGenerator
 *
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/GenerateDICOMDIR.exe path output_filename
 */
using System;
using gdcm;

public class GenerateDICOMDIR
{
    public static int Main(string[] args)
    {
        string directory = args[0];
        string outfilename = args[1];

        Directory d = new Directory();
        uint nfiles = d.Load( directory, true );
        if(nfiles == 0) return 1;
        //System.Console.WriteLine( "Files:\n" + d.toString() );

        // Implement fast path ?
        // Scanner s = new Scanner();

        string descriptor = "My_Descriptor";
        FilenamesType filenames = d.GetFilenames();

        gdcm.DICOMDIRGenerator gen = new DICOMDIRGenerator();
        gen.SetFilenames( filenames );
        gen.SetDescriptor( descriptor );
        if( !gen.Generate() )
        {
            return 1;
        }

        gdcm.FileMetaInformation.
            SetSourceApplicationEntityTitle( "GenerateDICOMDIR" );
        gdcm.Writer writer = new Writer();
        writer.SetFile( gen.GetFile() );
        writer.SetFileName( outfilename );
        if( !writer.Write() )
        {
            return 1;
        }

        return 0;
    }
}

```

## 29.73 GenerateRTSTRUCT.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre

```

All rights reserved.  
See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even  
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR  
PURPOSE. See the above copyright notice for more information.

```

=====*/
#include "vtkGDCMPolyDataWriter.h"
#include "vtkGDCMPolyDataReader.h"
#include "vtkPolyData.h"
#include "vtkPolyDataReader.h"
#include "vtkMedicalImageProperties.h"
#include "vtkRTStructSetProperties.h"
#include "vtkStringArray.h"
#include "vtkAppendPolyData.h"
#include "vtkPolyDataWriter.h"
#include "vtkPolyDataMapper.h"
#include "vtkPolyDataMapper2D.h"
#include "vtkActor2D.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkMedicalImageProperties.h"
#include "vtkRenderWindow.h"
#include "vtkRenderer.h"
#include "vtkCamera.h"
#include "vtkProperty.h"
#include "vtkProperty2D.h"
#include "vtkImageData.h"

#include <algorithm> //for std::find

#include "gdcmDirectoryHelper.h"

using namespace gdcm;

//view each organ independently of the others, to make sure that
//organ names correspond to actual segmentations.
void ShowOrgan(vtkPolyData* inData)
{
    // Now we'll look at it.
    vtkPolyDataMapper *cubeMapper = vtkPolyDataMapper::New();
    #if (VTK_MAJOR_VERSION >= 6)
        cubeMapper->SetInputData( inData );
    #else
        cubeMapper->SetInput( inData );
    #endif
    cubeMapper->SetScalarRange(0,7);
    vtkActor *cubeActor = vtkActor::New();
    cubeActor->SetMapper(cubeMapper);
    vtkProperty * property = cubeActor->GetProperty();
    property->SetRepresentationToWireframe();

    vtkRenderer *renderer = vtkRenderer::New();
    vtkRenderWindow *renWin = vtkRenderWindow::New();
    renWin->AddRenderer(renderer);

    vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
    iren->SetRenderWindow(renWin);

    renderer->AddActor(cubeActor);
    renderer->ResetCamera();
    renderer->SetBackground(1,1,1);

    renWin->SetSize(300,300);

    renWin->Render();
    iren->Start();

    cubeMapper->Delete();
    cubeActor->Delete();
    renderer->Delete();
    renWin->Delete();
    iren->Delete();
}

/*
 * Full application which ... RTSTRUCT
 */
int main(int argc, char *argv[])
{
    if( argc < 2 )

```

```

    {
        std::cerr << argv[0] << " directory-with-rtstruct-and-ct-images\n";
        return 1;
    }
    std::string theDirName(argv[1]);
    Directory::FileNamesType theRTSeries =
        DirectoryHelper::GetRTStructSeriesUIDs(theDirName);

    gdcmm::Directory theDir;
    theDir.Load(argv[1]);

    if (theRTSeries.empty())
    {
        std::cerr << "No RTStructs found for the test, ending." << std::endl;
        return 1;
    }

    for (size_t q = 0; q < theRTSeries.size(); q++)
    {
        Directory::FileNamesType theRTNames =
            DirectoryHelper::GetFileNamesFromSeriesUIDs(theDirName,
                theRTSeries[q]);

        if (theRTNames.empty()) {
            std::cerr << "Unable to load RT Series " << theRTSeries[q] << ", continuing. " << std::endl;
            continue;
        }

        vtkGDCMPolyDataReader * reader =
            vtkGDCMPolyDataReader::New();
        reader->SetFileName( theRTNames[0].c_str() );
        reader->Update();

        //std::cout << reader->GetMedicalImageProperties()->GetStudyDate() << std::endl;

        vtkGDCMPolyDataWriter * writer =
            vtkGDCMPolyDataWriter::New();
        int numMasks = reader->GetNumberOfOutputPorts() + 1; //add a blank one in
        writer->SetNumberOfInputPorts( numMasks );
        std::string thePotentialName = theDirName + "/" + "GDCMTestRTStruct." + theRTSeries[q] + ".dcm";
        gdcmm::Directory::FileNamesType theFileNames = theDir.
            GetFileNames();
        //keep renaming the output until we get something that doesn't overwrite what was there already
        int count = 0;
        while (std::find(theFileNames.begin(), theFileNames.end(), thePotentialName) != theFileNames.end())
        {
            char buff[255];
            sprintf(buff, "%d", count);
            thePotentialName = theDirName + "/" + "GDCMTestRTStruct." + buff + "." + theRTSeries[q] + ".dcm";
        }
        writer->SetFileName( thePotentialName.c_str() );
        writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
        //this line is cheating, we won't have the same stuff, and may not have a struct
        //to start with.
        //have to go back to the original data to reconstruct the RTStructureSetProperties
        //writer->SetRTStructureSetProperties( reader->GetRTStructureSetProperties() );
        //writer->Write();

        //loop through the outputs in order to write them out as if they had been created and appended
        vtkStringArray* roiNames = vtkStringArray::New();
        vtkStringArray* roiAlgorithms = vtkStringArray::New();
        vtkStringArray* roiTypes = vtkStringArray::New();
        roiNames->SetNumberOfValues(numMasks);
        roiAlgorithms->SetNumberOfValues(numMasks);
        roiTypes->SetNumberOfValues(numMasks);
        vtkAppendPolyData* append = vtkAppendPolyData::New();

        //ok, now we'll add a blank organ
        //the blank organ is to test to ensure that blank organs work; there have been crash reports
        //this code is added at the beginning to ensure that the blank organs are read
        //and preserved as individual organs.
        vtkPolyData* blank = vtkPolyData::New();
        #if (VTK_MAJOR_VERSION >= 6)
            writer->SetInputData(0, blank);
        #else
            writer->SetInput(0, blank);
        #endif
        roiNames->InsertValue(0, "blank");
        roiAlgorithms->InsertValue(0, "blank");
        roiTypes->InsertValue(0, "ORGAN");
    }

```

```

//note the offsets used to place the blank rtstruct at the beginning of the newly generated RT.
//the idea is to run the program twice; first to generate an rtstruct with a blank mask (making
//sure that that functionality works), and then a second time to make sure that everything is
//being read properly. Multiple organs with the same name could cause some strangenesses.
for (int i = 1; i < numMasks; ++i)
{
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputConnection(i, reader->GetOutputPort(i-1));
        append->AddInputConnection(reader->GetOutputPort(i-1));
    #else
        writer->SetInput(i, reader->GetOutput(i-1));
        append->AddInput(reader->GetOutput(i-1));
    #endif

    std::string theString = reader->GetRTStructSetProperties()->GetStructureSetROIName(i-1);
    roiNames->InsertValue(i, theString);
    theString = reader->GetRTStructSetProperties()->GetStructureSetROIGenerationAlgorithm(i-1);
    roiAlgorithms->InsertValue(i, theString);
    theString = reader->GetRTStructSetProperties()->GetStructureSetRTROIInterpretedType(i-1);
    roiTypes->InsertValue(i, theString);

    ShowOrgan(reader->GetOutput(i-1));
}

vtkRTStructSetProperties* theProperties =
    vtkRTStructSetProperties::New();
writer->SetRTStructSetProperties(theProperties);
writer->InitializeRTStructSet(theDirName,
    reader->GetRTStructSetProperties()->GetStructureSetLabel(),
    reader->GetRTStructSetProperties()->GetStructureSetName(),
    roiNames, roiAlgorithms, roiTypes);

writer->SetRTStructSetProperties(theProperties);
writer->Write();

// print reader output:
reader->Print( std::cout );
// print first output:
reader->GetOutput()->Print( std::cout );

reader->Delete();
append->Delete();
roiNames->Delete();
roiTypes->Delete();
theProperties->Delete();
roiAlgorithms->Delete();
blank->Delete();

writer->Delete();
}
return 0;
}

```

## 29.74 GenerateStandardSOPClasses.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
*/

#include "gdcmDefs.h"
#include "gdcmUIDs.h"
#include "gdcmGlobal.h"
#include "gdcmMediaStorage.h"
#include "gdcmSOPClassUIDToIOD.h"

int main(int , char *[])

```

```

{
    using gdcm::MediaStorage;
    gdcm::Global& g = gdcm::Global::GetInstance();
    if( !g.LoadResourcesFiles() )
    {
        std::cerr << "Could not LoadResourcesFiles" << std::endl;
        return 1;
    }

    const gdcm::Defs &defs = g.GetDefs();

    int ret = 0;

    //std::cout << "Table B.5-1 STANDARD SOP CLASSES" << std::endl;
    std::cout << "SOP Class Name,SOP Class UID,IOD Specification (defined in PS 3.3)" << std::endl;

    gdcm::MediaStorage::MSType mst;
    for ( mst = gdcm::MediaStorage::MediaStorageDirectoryStorage
        ; mst < gdcm::MediaStorage::MS_END;
        mst = (gdcm::MediaStorage::MSType)(mst + 1) )
    {
        const char *iod = defs.GetIODNameFromMediaStorage(mst);
        gdcm::UIDs uid;
        uid.SetFromUID( gdcm::MediaStorage::GetMSString(mst) /*
            mst.GetString() */ );
        if( iod )
        {
            const char *iod_ref = gdcm::SOPClassUIDToIOD::GetIOD(uid);
            if( iod_ref )
            {
                std::string iod_ref_str = iod_ref;
                //iod_ref_str += " IOD Modules";
                //if( iod_ref_str != iod )
                {
                    //std::cout << "UID: " << uid << " ";
                    std::cout << "' ' << uid.GetName() << "' ' << ", " << "' ' << uid.
                        GetString() << "' ' << ", " << "' ' << iod << "' ' << std::endl;
                    //std::cout << "Incompatible IODs: [" << iod << "]" versus ref= [" << iod_ref_str << "]" <<
                        std::endl;
                    ++ret;
                }
            }
        }
    }

    return 0;
}

```

## 29.75 GenFakelIdentifyFile.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmGlobal.h"
#include "gdcmDummyValueGenerator.h"
#include "gdcmMediaStorage.h"
#include "gdcmWriter.h"
#include "gdcmItem.h"
#include "gdcmImageReader.h"
#include "gdcmSequenceOfItems.h"
#include "gdcmAttribute.h"
#include "gdcmFile.h"
#include "gdcmTag.h"
#include "gdcmDict.h"

```



```

#include "gdcDictEntry.h"
#include "gdcDicts.h"
#include "gdcTransferSyntax.h"
#include "gdcUIDGenerator.h"
#include "gdcAnonymizer.h"

#include <cstdlib>
#include <cstring>

gdc::DataElement CreateFakeElement(gdc::Tag const &tag, bool toremove)
{
    static const gdc::Global &g = gdc::Global::GetInstance();
    static const gdc::Dicts &dicts = g.GetDicts();
    static const gdc::Dict &pubdict = dicts.GetPublicDict();
    static size_t countglobal = 0;
    static std::vector<gdc::Tag> balcptags =
        gdc::Anonymizer::GetBasicApplicationLevelConfidentialityProfileAttributes
        ();
    size_t count = countglobal % balcptags.size();

    const gdc::DictEntry &dictentry = pubdict.GetDictEntry(tag);

    gdc::DataElement de;
    de.SetTag( tag );
    using gdc::VR;
    const VR &vr = dictentry.GetVR();
    //if( vr != VR::INVALID )
    if( vr.IsDual() )
    {
        if( vr == VR::US_SS )
        {
            de.SetVR( VR::US );
        }
        else if( vr == VR::US_SS_OW )
        {
            de.SetVR( VR::OW );
        }
        else if( vr == VR::OB_OW )
        {
            de.SetVR( VR::OB );
        }
    }
    else
    {
        de.SetVR( vr );
    }
    const char str[] = "BasicApplicationLevelConfidentialityProfileAttributes";
    const char safe[] = "This is safe to keep";
    if( de.GetVR() != VR::SQ )
    {
        if( toremove )
            de.SetByteValue( str, (uint32_t)strlen(str) );
        else
            de.SetByteValue( safe, (uint32_t)strlen(safe) );
    }
    else
    {
        // Create an item
        gdc::Item it;
        it.SetVLToUndefined();
        gdc::DataSet &nds = it.GetNestedDataSet();
        // Insert sequence into data set
        assert(de.GetVR() == gdc::VR::SQ );
        gdc::SmartPointer<gdc::SequenceOfItems> sq = new
            gdc::SequenceOfItems();
        sq->SetLengthToUndefined();
        de.SetValue(*sq);
        de.SetVLToUndefined();
        //ds.Insert( de);

        if( !toremove )
        {
            nds.Insert( CreateFakeElement( balcptags[count], true ) );
            countglobal++;
        }
        else
        {
            gdc::Attribute<0x0008,0x0000> at1 = { 0 }; // This element has no
                reason to be 'anonymized'...
            nds.Insert( at1.GetAsDataElement() );
            gdc::Attribute<0x000a,0x0000> at2 = { 0 };
        }
    }
}

```

```

        nds.Insert( at2.GetAsDataElement() );
    }
    sq->AddItem(it);
}
return de;
}

/*
*/
int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " output.dcm" << std::endl;
        return 1;
    }
    using gdcm::Tag;
    using gdcm::VR;
    const char *outfilename = argv[1];

    std::vector<gdcm::Tag> balcptags =
        gdcm::Anonymizer::GetBasicApplicationLevelConfidentialityProfileAttributes
            ();

    gdcm::Writer w;
    gdcm::File &f = w.GetFile();
    gdcm::DataSet &ds = f.GetDataSet();

    // Add attribute that need to be anonymized:
    std::vector<gdcm::Tag>::const_iterator it = balcptags.begin();
    for(; it != balcptags.end(); ++it)
    {
        ds.Insert( CreateFakeElement( *it, true ) );
    }

    // Add attribute that do NOT need to be anonymized:
    static const gdcm::Global &g = gdcm::Global::GetInstance();
    static const gdcm::Dicts &dicts = g.GetDicts();
    static const gdcm::Dict &pubdict = dicts.GetPublicDict();

    using gdcm::Dict;
    Dict::ConstIterator dictit = pubdict.Begin();
    for(; dictit != pubdict.End(); ++dictit)
    {
        const gdcm::Tag &dicttag = dictit->first;
        if( dicttag == Tag(0x6e65,0x6146) ) break;
        //const gdcm::DictEntry &dictentry = dictit->second;
        ds.Insert( CreateFakeElement( dicttag, false ) );
    }
    ds.Remove( gdcm::Tag(0x400,0x500) );
    ds.Remove( gdcm::Tag(0x12,0x62) );
    ds.Remove( gdcm::Tag(0x12,0x63) );

    // Make sure to override any UID stuff
    gdcm::UIDGenerator uid;
    gdcm::DataElement de( Tag(0x8,0x18) ); // SOP Instance UID
    de.SetVR( VR::UI );
    const char *u = uid.Generate();
    de.SetByteValue( u, (uint32_t)strlen(u) );
    //ds.Insert( de );
    ds.Replace( de );

    de.SetTag( Tag(0x8,0x16) ); // SOP Class UID
    de.SetVR( VR::UI );
    gdcm::MediaStorage ms( gdcm::MediaStorage::RawDataStorage
        );
    de.SetByteValue( ms.GetString(), (uint32_t)strlen(ms.GetString()) );
    ds.Replace( de ); // replace !

    gdcm::FileMetaInformation &fmi = f.GetHeader();
    //fmi.SetDataSetTransferSyntax( gdcm::TransferSyntax::ImplicitVRLittleEndian );
    fmi.SetDataSetTransferSyntax(
        gdcm::TransferSyntax::ExplicitVRLittleEndian );

    w.SetCheckFileMetaInformation( true );
    w.SetFileName( outfile );
    if (!w.Write() )
    {
        return 1;
    }
}

```

```

    return 0;
}

```

## 29.76 GenFakelImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmImage.h"
#include "gdcmImageWriter.h"
#include "gdcmFileDerivation.h"
#include "gdcmUIDGenerator.h"
//#include "gdcmImageChangePhotometricInterpretation.h"

/*
 * This example shows two things:
 * 1. How to create an image ex-nihilo
 * 2. How to use the gdcm.FileDerivation filter. This filter is meant to create "DERIVED" image
 * object. FileDerivation has a simple API where you can reference *all* the input image that have been
 * used to generate the image. The API also allows user to specify the purpose of reference (see CID 7202,
 * PS 3.16 - 2008), and the image derivation type (CID 7203, PS 3.16 - 2008).
 */
int main(int, char *[])
{
    // Step 1: Fake Image
    gdcm::SmartPointer<gdcm::Image> im = new
        gdcm::Image;

    char * buffer = new char[ 256 * 256 * 3];
    char * p = buffer;
    int b = 128;
    //int ybr[3];
    int ybr2[3];
    //int rgb[3];

    for(int r = 0; r < 256; ++r)
        for(int g = 0; g < 256; ++g)
            //for(int b = 0; b < 256; ++b)
            {
                //rgb[0] = r;
                //rgb[1] = g;
                //rgb[1] = 128;
                //rgb[2] = b;
                //ybr[0] = r;
                //ybr[1] = g;
                //ybr[1] = 128;
                //ybr[2] = b;

                ybr2[0] = r;
                ybr2[1] = g;
                ybr2[1] = 128;
                ybr2[2] = b;
                //gdcm::ImageChangePhotometricInterpretation::YBR2RGB(rgb, ybr);
                //gdcm::ImageChangePhotometricInterpretation::RGB2YBR(ybr2, rgb);
                *p++ = (char)ybr2[0];
                *p++ = (char)ybr2[1];
                *p++ = (char)ybr2[2];
            }

    im->SetNumberOfDimensions( 2 );
    im->SetDimension(0, 256 );
    im->SetDimension(1, 256 );

    im->GetPixelFormat().SetSamplesPerPixel(3);
    //im->SetPhotometricInterpretation( gdcm::PhotometricInterpretation::RGB );
    im->SetPhotometricInterpretation(
        gdcm::PhotometricInterpretation::YBR_FULL );

```

```

unsigned long l = im->GetBufferLength();
if( l != 256 * 256 * 3 )
{
    return 1;
}
gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
pixeldata.SetByteValue( buffer, (uint32_t)l );
delete[] buffer;
im->SetDataElement( pixeldata );

gdcm::UIDGenerator uid; // helper for uid generation

gdcm::SmartPointer<gdcm::File> file = new
    gdcm::File; // empty file

// Step 2: DERIVED object
gdcm::FileDerivation fd;
// For the pupose of this exercise we will pretend that this image is referencing
// two source image (we need to generate fake UID for that).
const char ReferencedSOPClassUID[] = "1.2.840.10008.5.1.4.1.1.7"; // Secondary Capture
fd.AddReference( ReferencedSOPClassUID, uid.Generate() );
fd.AddReference( ReferencedSOPClassUID, uid.Generate() );

// Again for the purpose of the exercise we will pretend that the image is a
// multiplanar reformat (MPR):
// CID 7202 Source Image Purposes of Reference
// {"DCM",121322,"Source image for image processing operation"},
fd.SetPurposeOfReferenceCodeSequenceCodeValue( 121322 );
// CID 7203 Image Derivation
// { "DCM",113072,"Multiplanar reformatting" },
fd.SetDerivationCodeSequenceCodeValue( 113072 );
fd.SetFile( *file );
// If all Code Value are ok the filter will execute properly
if( !fd.Derive() )
{
    std::cerr << "Sorry could not derive using input info" << std::endl;
    return 1;
}

// We pass both :
// 1. the fake generated image
// 2. the 'DERIVED' dataset object
// to the writer.
gdcm::ImageWriter w;
w.SetImage( *im );
w.SetFile( fd.GetFile() );

// Set the filename:
w.SetFileName( "ybr2.dcm" );
if( !w.Write() )
{
    return 1;
}

return 0;
}

```

## 29.77 GenLongSeqs.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmItem.h"
#include "gdcmImageReader.h"

```

```

#include "gdcmSequenceOfItems.h"
#include "gdcmFile.h"
#include "gdcmTag.h"

/*
 * This example is used to generate the file:
 *
 *
 * There is a flaw in the DICOM design where it is assumed that Sequence can be
 * either represented as undefined length or defined length. This should work
 * in most case, but the undefined length is a little more general and can
 * store sequence of items that a defined length cannot.
 * We need to make sure that we can store numerous Item in a SQ
 *
 * Warning: do not try to compute the group length elements !
 * Warning: You may need a 64bits machine for this example to work.
 */
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();

    // Create a Sequence
    gdcm::SmartPointer<gdcm::SequenceOfItems> sq = new
        gdcm::SequenceOfItems();
    sq->SetLengthToUndefined();

    const char owner_str[] = "GDCM CONFORMANCE TESTS";
    gdcm::DataElement owner( gdcm::Tag(0x4d4d, 0x10) );
    owner.SetByteValue( owner_str, (uint32_t)strlen(owner_str));
    owner.SetVR( gdcm::VR::LO );

    size_t nitems = 1000;
    nitems += std::numeric_limits<uint32_t>::max();
    for(unsigned int idx = 0; idx < nitems; ++idx)
    {
        // Create a dataelement
        //gdcm::DataElement de( gdcm::Tag(0x4d4d, 0x1002) );
        //de.SetByteValue(ptr, ptr_len);
        //de.SetVR( gdcm::VR::OB );

        // Create an item
        gdcm::Item it;
        it.SetVLToUndefined();
        //gdcm::DataSet &nds = it.GetNestedDataSet();
        //nds.Insert(owner);
        //nds.Insert(de);

        sq->AddItem(it);
    }

    // Insert sequence into data set
    gdcm::DataElement des( gdcm::Tag(0x4d4d, 0x1001) );
    des.SetVR(gdcm::VR::SQ);
    des.SetValue(*sq);
    des.SetVLToUndefined();

    ds.Insert( owner );
    ds.Insert( des );

    gdcm::Writer w;
    w.SetFile( file );
    //w.SetCheckFileMetaInformation( true );
    w.SetFileName( outfile );
    if( !w.Write() )
    {
        return 1;
    }
}

```

```

    }
    return 0;
}

```

## 29.78 GenSeqs.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
#include "gdcmlReader.h"
#include "gdcmlWriter.h"
#include "gdcmlItem.h"
#include "gdcmlImageReader.h"
#include "gdcmlSequenceOfItems.h"
#include "gdcmlFile.h"
#include "gdcmlTag.h"

/*
 * This example is used to generate the file:
 *
 * gdcmlConformanceTests/SequenceWithUndefinedLengthNotConvertibleToDefinedLength.dcm
 *
 * There is a flaw in the DICOM design where it is assumed that Sequence can be
 * either represented as undefined length or defined length. This should work
 * in most cases, but the undefined length is a little more general and can
 * store sequence of items that a defined length cannot.
 * Deflated syntax was used in this case since this synthetic example can be
 * nicely compressed using this transfer syntax.
 *
 * Warning: do not try to compute the group length elements !
 * Warning: You may need a 64bits machine for this example to work.
 */
int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcml::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcml::File &file = reader.GetFile();
    gdcml::DataSet &ds = file.GetDataSet();

    //const unsigned int nitems = 1000;
    const unsigned int ptr_len = 42; /*94967296 / nitems; */
    //assert( ptr_len == 42949672 );
    char *ptr = new char[ptr_len];
    memset(ptr,0,ptr_len);

    // Create a Sequence
    gdcml::SmartPointer<gdcml::SequenceOfItems> sq = new
        gdcml::SequenceOfItems();
    sq->SetLengthToUndefined();

    const char owner_str[] = "GDCM CONFORMANCE TESTS";
    gdcml::DataElement owner( gdcml::Tag(0x4d4d, 0x10) );
    owner.SetByteValue(owner_str, (uint32_t)strlen(owner_str));
    owner.SetVR( gdcml::VR::LO );

```

```

for(unsigned int idx = 0; idx < 10/* nitems*/; ++idx)
{
    // Create a dataelement
    gdcm::DataElement de( gdcm::Tag(0x4d4d, 0x1002) );
    de.SetByteValue(ptr, ptr_len);
    de.SetVR( gdcm::VR::OB );

    // Create an item
    gdcm::Item it;
    it.SetVLToUndefined();
    gdcm::DataSet &nds = it.GetNestedDataSet();
    nds.Insert(owner);
    nds.Insert(de);

    sq->AddItem(it);
}

// Insert sequence into data set
gdcm::DataElement des( gdcm::Tag(0x4d4d,0x1001) );
des.SetVR(gdcm::VR::SQ);
des.SetValue(*sq);
des.SetVLToUndefined();

ds.Insert(owner);
ds.Insert(des);

gdcm::Writer w;
w.SetFile( file );
//w.SetCheckFileMetaInformation( true );
w.SetFileName( outfilename );
if (!w.Write() )
{
    return 1;
}

return 0;
}

```

## 29.79 GetArray.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/GetArray.exe gdcmData/012345.002.050.dcm
 */
using System;
using gdcm;

public class GetArray
{
    public static int Main(string[] args)
    {
        {
            string file1 = args[0];
            ImageReader reader = new ImageReader();
            reader.SetFileName( file1 );
            bool ret = reader.Read();
            if( !ret )
            {
                return 1;
            }

            Image image = reader.GetImage();

```

```

PixelFormat pixeltype = image.GetPixelFormat();

if( image.GetNumberOfDimensions() != 2 )
{
    // For the purpose of the test, exit early on
    return 1;
}
uint dimx = image.GetDimension(0);
uint dimy = image.GetDimension(1);
uint npixels = dimx * dimy;
//LookupTable lut = image.GetLUT();
//uint r1 = lut.GetLUTLength( LookupTable.LookupTableType.RED );
//byte[] rbuf = new byte[ r1 ];
//uint r12 = lut.GetLUT( LookupTable.LookupTableType.RED, rbuf );
//assert r1 == r12;

//byte[] str1 = new byte[ image.GetBufferLength()];
//image.GetBuffer( str1 );
if( pixeltype.GetScalarType() == PixelFormat.ScalarType.UINT8 )
{
    System.Console.WriteLine( "Processing UINT8 image type" );
    byte[] str1 = new byte[ npixels ];
    image.GetArray( str1 );
}
else if( pixeltype.GetScalarType() == PixelFormat.ScalarType.INT16 )
{
    System.Console.WriteLine( "Processing INT16 image type" );
    short[] str1 = new short[ npixels ];
    image.GetArray( str1 );
}
else if( pixeltype.GetScalarType() == PixelFormat.ScalarType.UINT16 )
{
    System.Console.WriteLine( "Processing UINT16 image type" );
    ushort[] str1 = new ushort[ npixels ];
    image.GetArray( str1 );
}
else
{
    //System.Console.WriteLine( "Default (unhandled pixel format): " + pixeltype.ToString() );
    System.Console.WriteLine( "Default (unhandled pixel format): " + pixeltype.GetScalarTypeAsString() );
    // Get bytes
    byte[] str1 = new byte[ image.GetBufferLength()];
    image.GetBuffer( str1 );
}

return 0;
}
}

```

## 29.80 GetJPEGSamplePrecision.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
* This example is a little helper to detect the famous SIEMENS JPEG lossless compressed image
* where DICOM is declared as:
*
* (0028,0100) US 16 # 2,1 Bits Allocated
* (0028,0101) US 12 # 2,1 Bits Stored
* (0028,0102) US 11 # 2,1 High Bit
* (0028,0103) US 0 # 2,1 Pixel Representation
*
* But where JPEG is:
*
* JPEG_SOF_Parameters:

```



```

*           SamplePrecision = 16
*           nLines = 192
*           nSamplesPerLine = 192
*           nComponentsInFrame = 1
*           component 0
*               ComponentIdentifier = 1
*               HorizontalSamplingFactor = 1
*               VerticalSamplingFactor = 1
*               QuantizationTableDestinationSelector = 0
*
* This case is valid. One simply has to use the 16bits jpeg decoder to decode the 12bits stored image.
* This used to be an issue in GDCM 1.2.x (fixed in GDCM 1.2.5)
*
* The main return 0 (no error) when the file read is actually a potential problem. At the end of the main
* function, the jpeg stream is stored in the filename specified as second argument
*/

#include "gdcmImageReader.h"
#include "gdcmSequenceOfFragments.h"
#include "gdcmJPEGCodec.h"

#include <iostream>
#include <fstream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.jpg" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    gdcm::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }

    // The output of gdcm::Reader is a gdcm::File
    const gdcm::File &file = reader.GetFile();
    const gdcm::Image &image = reader.GetImage();

    const gdcm::TransferSyntax &ts = file.GetHeader().
        GetDataSetTransferSyntax();

    if( ts != gdcm::TransferSyntax::JPEGLosslessProcess14 && ts !=
        gdcm::TransferSyntax::JPEGLosslessProcess14_1 )
    {
        std::cerr << "Input is not a lossless JPEG" << std::endl;
        return 1;
    }

    // the dataset is the the set of element we are interested in:
    const gdcm::DataSet &ds = file.GetDataSet();

    const gdcm::Tag rawTag(0x7fe0, 0x0010); // Default to Pixel Data
    const gdcm::DataElement& pdde = ds.GetDataElement( rawTag );
    const gdcm::SequenceOfFragments *sf = pdde.
        GetSequenceOfFragments();
    if( sf )
    {
        std::ofstream output(outfilename, std::ios::binary);
        sf->WriteBuffer(output);
    }
    else
    {
        std::cerr << "Error" << std::endl;
        return 1;
    }

    gdcm::JPEGCodec jpeg;
    std::ifstream is(outfilename, std::ios::binary);
    gdcm::PixelFormat pf ( gdcm::PixelFormat::UINT8 ); // let's
        pretend it's a 8bits jpeg
    jpeg.SetPixelFormat( pf );
    gdcm::TransferSyntax ts_jpg;

```

```

bool b = jpeg.GetHeaderInfo( is, ts_jpg );
if( !b )
{
    return 1;
}

//jpeg.Print( std::cout );
if( jpeg.GetPixelFormat().GetBitsAllocated() != image.
    GetPixelFormat().GetBitsAllocated()
|| jpeg.GetPixelFormat().GetBitsStored() != image.
    GetPixelFormat().GetBitsStored() )
{
    std::cerr << "There is a mismatch in between DICOM declared Pixel Format and Sample Precision used in
    the JPEG stream" << std::endl;
    return 0;
}

std::cout << jpeg.GetPixelFormat() << std::endl;
std::cout << image.GetPixelFormat() << std::endl;

return 1;
}

```

## 29.81 GetPortionCSAHeader.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python GetPortionCSAHeader.py input.dcm
19
20 Footnote:
21 SIEMENS is not publishing any information on the CSA header. So any info extracted
22 is at your own risk.
23 """
24
25 import sys
26 import gdcm
27
28 if __name__ == "__main__":
29
30     file = sys.argv[1]
31
32     r = gdcm.Reader()
33     r.SetFileName( file )
34     if not r.Read():
35         sys.exit(1)
36
37     ds = r.GetFile().GetDataSet()
38     csa_t1 = gdcm.CSAHeader()
39     csa_t2 = gdcm.CSAHeader()
40     #print csa
41     t1 = csa_t1.GetCSAImageHeaderInfoTag();
42     print t1
43     t2 = csa_t2.GetCSASeriesHeaderInfoTag();
44     print t2
45     # Let's do it for t1:
46     if ds.FindDataElement( t1 ):
47         csa_t1.LoadFromDataElement( ds.GetDataElement( t1 ) )
48         print csa_t1
49
50 # Now let's pretend we are only interested in B_value and DiffusionGradientDirection entries:
51 bvalues = csa_t1.GetCSAElementByName( "B_value" ) # WARNING: it is case sensitive !
52 print bvalues

```

```

53
54   diffgraddir = csa_t1.GetCSAElementByName( "DiffusionGradientDirection" ) # WARNING: it is case sensitive
55   !
56   print diffgraddir
57   # repeat for t2 if you like it:
58   if ds.FindDataElement( t2 ):
59       csa_t2.LoadFromDataElement( ds.GetDataElement( t2 ) )
60       # print csa_t2
61
62   gdt = csa_t2.GetCSAElementByName( "GradientDelayTime" )
63   print gdt
64
65   bv = gdt.GetByteValue();
66   #print bv
67   str = bv.GetPointer()
68   print str.split("\\")

```

## 29.82 GetSequenceUltrasound.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmAttribute.h"

bool Region ( char* nomefile, unsigned int* X_min, unsigned int* Y_min, unsigned int* X_max, unsigned int*
    Y_max );

int main(int argc, char* argv[] )
{
    // Controllo del numero di argomenti introdotti da riga di comando
    if( argc < 2 )
    {
        std::cerr << "Usage: " << std::endl;
        std::cerr << argv[0] << " inputImageFile " << std::endl;
        return EXIT_FAILURE;
    }

    unsigned int x_min = 1;
    unsigned int y_min = 1;
    unsigned int x_max = 1;
    unsigned int y_max = 1;

    if( Region ( argv[1], &x_min, &y_min, &x_max, &y_max ) )
    {
        std::cout << "x_min = " << x_min << std::endl;
        std::cout << "y_min = " << y_min << std::endl;
        std::cout << "x_max = " << x_max << std::endl;
        std::cout << "y_max = " << y_max << std::endl;
    }

    else
    {
        std::cout << "no\n";
    }
}

bool Region ( char* nomefile, unsigned int* X_min, unsigned int* Y_min, unsigned int* X_max, unsigned int*
    Y_max )
{
    gdcm::Reader reader;
    reader.SetFileName( nomefile );
    if( !reader.Read() )
    {

```

```

    std::cerr << "Could not read: " << nomefile << std::endl;
    return false;
}

gdcmm::File &file = reader.GetFile();
gdcmm::DataSet &ds = file.GetDataSet();

gdcmm::Tag tsqr(0x0018,0x6011);
if( !ds.FindDataElement( tsqr ) )
{
    return false;
}

const gdcmm::DataElement &sqr= ds.GetDataElement( tsqr );
//std::cout << sqr << std::endl;
const gdcmm::SequenceOfItems *sqi = sqr.GetValueAssQ();
if( !sqi || !sqi->GetNumberOfItems() )
{
    return false;
}
//std::cout << sqi << std::endl;

const gdcmm::Item &item = sqi->GetItem(1);
//std::cout << item << std::endl;
const gdcmm::DataSet& nestedds = item.GetNestedDataSet();
//std::cout << nestedds << std::endl;

gdcmm::Tag tX0(0x0018,0x6018);
gdcmm::Tag tY0(0x0018,0x601a);
gdcmm::Tag tX1(0x0018,0x601c);
gdcmm::Tag tY1(0x0018,0x601e);

if( (!nestedds.FindDataElement( tX0 ))||(!nestedds.
    FindDataElement( tY0 ))||(!nestedds.FindDataElement( tX1 ))||(!nestedds.
    FindDataElement( tY1 )) )
{
    return false;
}

const gdcmm::DataElement& deX0 = nestedds.GetDataElement( tX0 );
const gdcmm::DataElement& deY0 = nestedds.GetDataElement( tY0 );
const gdcmm::DataElement& deX1 = nestedds.GetDataElement( tX1 );
const gdcmm::DataElement& deY1 = nestedds.GetDataElement( tY1 );
//std::cout << deX0 << std::endl << deY0 << std::endl << deX1 << std::endl << deY1 << std::endl;

//const gdcmm::ByteValue *bvX0 = deX0.GetByteValue();
//const gdcmm::ByteValue *bvY0 = deY0.GetByteValue();
//const gdcmm::ByteValue *bvX1 = deX1.GetByteValue();
//const gdcmm::ByteValue *bvY1 = deY1.GetByteValue();
//std::cout << bvX0 << std::endl << bvY0 << std::endl << bvX1 << std::endl << bvY1 << std::endl;

gdcmm::Attribute<0x0018,0x6018> atX0;
gdcmm::Attribute<0x0018,0x601a> atY0;
gdcmm::Attribute<0x0018,0x601c> atX1;
gdcmm::Attribute<0x0018,0x601e> atY1;
atX0.SetFromDataElement( deX0 );
atY0.SetFromDataElement( deY0 );
atX1.SetFromDataElement( deX1 );
atY1.SetFromDataElement( deY1 );
uint32_t X0 = atX0.GetValue();
uint32_t Y0 = atY0.GetValue();
uint32_t X1 = atX1.GetValue();
uint32_t Y1 = atY1.GetValue();
std::cout << X0 << std::endl << Y0 << std::endl << X1 << std::endl << Y1 << std::endl;

*X_min = static_cast<unsigned int>(X0);
*Y_min = static_cast<unsigned int>(Y0);
*X_max = static_cast<unsigned int>(X1);
*Y_max = static_cast<unsigned int>(Y1);

//std::cout << "X_min = " << *X_min << std::endl;
//std::cout << "Y_min = " << *Y_min << std::endl;
//std::cout << "X_max = " << *X_max << std::endl;
//std::cout << "Y_max = " << *Y_max << std::endl;

return true;
}

```

## 29.83 GetSubSequenceData.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmImage.h"
#include "gdcmImageWriter.h"
#include "gdcmDataElement.h"
#include "gdcmPrivateTag.h"
#include "gdcmUIDGenerator.h"

#include <iostream>
#include <string>

#include <map>

/*
 * This example will extract the Movie from the private group of
 * GEMS_Ultrasound_MovieGroup_001 See Attribute
 * (7fe1,60,GEMS_Ultrasound_MovieGroup_001)
 *
 * The output file will be stored in 'outvid.dcm' as
 * MultiframeGrayscaleByteSecondaryCaptureImageStorage
 */
int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    using namespace gdcm;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    reader.Read();

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();
    const PrivateTag tseq(0x7fe1,0x1,"GEMS_Ultrasound_MovieGroup_001");

    if( !ds.FindDataElement( tseq ) ) return 1;
    const DataElement& seq = ds.GetDataElement( tseq );

    SmartPointer<SequenceOfItems> sqi = seq.
        GetValueAsSQ();
    assert( sqi->GetNumberOfItems() == 1 );
    Item &item = sqi->GetItem(1);
    DataSet &subds = item.GetNestedDataSet();

    const PrivateTag tseq1(0x7fe1,0x10,"GEMS_Ultrasound_MovieGroup_001");

    if( !subds.FindDataElement( tseq1 ) ) return 1;
    const DataElement& seq1 = subds.GetDataElement( tseq1 );

    SmartPointer<SequenceOfItems> sqi2 = seq1.
        GetValueAsSQ();
    //int n = sqi2->GetNumberOfItems();
    int index = 1;
    Item &item2 = sqi2->GetItem(index);
    DataSet &subds2 = item2.GetNestedDataSet();

    const PrivateTag tseq2(0x7fe1,0x20,"GEMS_Ultrasound_MovieGroup_001");

    if( !subds2.FindDataElement( tseq2 ) ) return 1;
    const DataElement& seq2 = subds2.GetDataElement( tseq2 );

    // std::cout << seq2 << std::endl;

    SmartPointer<SequenceOfItems> sqi3 = seq2.
        GetValueAsSQ();
    size_t ni3 = sqi3->GetNumberOfItems(); (void)ni3;
    assert( sqi3->GetNumberOfItems() >= 1 );
}

```

```

Item &item3 = sqi3->GetItem(1);
DataSet &subds3 = item3.GetNestedDataSet();

const PrivateTag tseq6(0x7fel,0x26,"GEMS_Ultrasound_MovieGroup_001");
if( !subds3.FindDataElement( tseq6 ) ) return 1;
const DataElement& seq6 = subds3.GetDataElement( tseq6 );
SmartPointer<SequenceOfItems> sqi6 = seq6.
    GetValueAsSQ();
size_t ni6= sqi6->GetNumberOfItems();
assert( sqi6->GetNumberOfItems() >= 1 );
const PrivateTag tseq7(0x7fel,0x86,"GEMS_Ultrasound_MovieGroup_001");
int dimx = 0, dimy = 0;
for( size_t i6 = 1; i6 <= ni6; ++i6 )
{
    Item &item6 = sqi6->GetItem(i6);
    DataSet &subds6 = item6.GetNestedDataSet();

    if( subds6.FindDataElement( tseq7 ) )
    {
        Element<VR::SL, VM::VM4> el;
        el.SetFromDataElement( subds6.GetDataElement( tseq7 ) );
        std::cout << "El= " << el.GetValue() << std::endl;
        dimx = el.GetValue(0);
        dimy = el.GetValue(1);
    }
}

const PrivateTag tseq3(0x7fel,0x36,"GEMS_Ultrasound_MovieGroup_001");
if( !subds3.FindDataElement( tseq3 ) ) return 1;
const DataElement& seq3 = subds3.GetDataElement( tseq3 );

//      std::cout << seq3 << std::endl;

SmartPointer<SequenceOfItems> sqi4 = seq3.
    GetValueAsSQ();
size_t ni4= sqi4->GetNumberOfItems();
assert( sqi4->GetNumberOfItems() >= 1 );
const PrivateTag tseq8(0x7fel,0x37,"GEMS_Ultrasound_MovieGroup_001");
const PrivateTag tseq4(0x7fel,0x43,"GEMS_Ultrasound_MovieGroup_001");
const PrivateTag tseq5(0x7fel,0x60,"GEMS_Ultrasound_MovieGroup_001");

std::vector<char> imbuffer;
int dimz = 0;
for( size_t i4 = 1; i4 <= ni4; ++i4 )
{
    Item &item4 = sqi4->GetItem(i4);
    DataSet &subds4 = item4.GetNestedDataSet();

    if( !subds4.FindDataElement( tseq8 ) ) return 1;
    const DataElement& de8 = subds4.GetDataElement( tseq8 );
    Element<VR::UL, VM::VM1> ldimz;
    ldimz.SetFromDataElement( de8 );
    dimz += ldimz.GetValue();
    if( !subds4.FindDataElement( tseq4 ) ) return 1;
    const DataElement& seq4 = subds4.GetDataElement( tseq4 );
    if( !subds4.FindDataElement( tseq5 ) ) return 1;
    const DataElement& seq5 = subds4.GetDataElement( tseq5 );

    //      std::cout << seq4 << std::endl;
    //      std::cout << seq5 << std::endl;

    const ByteValue *bv4 = seq4.GetByteValue();
    (void)bv4;
#ifdef 0
    {
        std::ofstream out( "/tmp/mo4", std::ios::binary );
        out.write( bv4->GetPointer(), bv4->GetLength());
        out.close();
    }
#endif
    const ByteValue *bv5 = seq5.GetByteValue();
#ifdef 0
    {
        std::ofstream out( "/tmp/mo5", std::ios::binary );
        out.write( bv5->GetPointer(), bv5->GetLength());
        out.close();
    }
#endif

    std::cout << bv5->GetLength() << std::endl;
    imbuffer.insert( imbuffer.begin(), bv5->GetPointer(), bv5->

```

```

        GetPointer() + bv5->GetLength() );
    }
    DataElement fakedata;
    fakedata.SetByteValue( &imbuffer[0], (uint32_t)imbuffer.size() );

    gdcm::SmartPointer<gdcm::Image> im = new
        gdcm::Image;
    im->SetNumberOfDimensions( 3 );

    im->SetDimension(0, dimx );
    im->SetDimension(1, dimy );
    im->SetDimension(2, dimz );
    size_t l1 = imbuffer.size();
    (void)l1;
    size_t l2 = im->GetBufferLength();
    (void)l2;
    assert( im->GetBufferLength() == imbuffer.size() );
    im->SetPhotometricInterpretation(
        gdcm::PhotometricInterpretation::MONOCHROME2 );

    im->SetDataElement( fakedata );

    gdcm::ImageWriter w;
    w.SetImage( *im );
    DataSet &dataset = w.GetFile().GetDataSet();

    gdcm::UIDGenerator uid;
    gdcm::DataElement de( Tag(0x8,0x18) ); // SOP Instance UID
    de.SetVR( VR::UI );
    const char *u = uid.Generate();
    de.SetByteValue( u, (uint32_t)strlen(u) );
    //ds.Insert( de );
    dataset.Replace( de );

    de.SetTag( Tag(0x8,0x16) ); // SOP Class UID
    de.SetVR( VR::UI );
    gdcm::MediaStorage ms(
        gdcm::MediaStorage::MultiframeGrayscaleByteSecondaryCaptureImageStorage
    );
    de.SetByteValue( ms.GetString(), (uint32_t)strlen(ms.
        GetString()));
    dataset.Replace( de ); // replace !

    w.SetFileName( "outvid.dcm" );
    if( !w.Write() )
    {
        return 1;
    }

    return 0;
}

```

## 29.84 headsq2dcm.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17 python headsq2dcm.py -D /path/to/VTKData
18 """
19
20 import vtk
21 import vtkgdcm
22 from vtk.util.misc import vtkGetDataRoot

```

```

23 VTK_DATA_ROOT = vtkGetDataRoot()
24
25 reader = vtk.vtkVolume16Reader()
26 reader.SetDataDimensions(64, 64)
27 reader.SetDataByteOrderToLittleEndian()
28 reader.SetFilePrefix(VTK_DATA_ROOT + "/Data/headsq/quarter")
29 reader.SetImageRange(1, 93)
30 reader.SetDataSpacing(3.2, 3.2, 1.5)
31
32 cast = vtk.vtkImageCast()
33 cast.SetInput( reader.GetOutput() )
34 cast.SetOutputScalarTypeToUnsignedChar()
35
36 # By default this is creating a Multiframe Grayscale Word Secondary Capture Image Storage
37 writer = vtkgdcmm.vtkGDCMImageWriter()
38 writer.SetFileName( "headsq.dcm" )
39 writer.SetInput( reader.GetOutput() )
40 # cast -> Multiframe Grayscale Byte Secondary Capture Image Storage
41 #writer.SetInput( cast.GetOutput() )
42 writer.SetFileDimensionality( 3 )
43 writer.Write()

```

## 29.85 HelloActiviz.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using vtkgdcmm;
using Kitware.VTK;
using System;
using System.Runtime.InteropServices;

/*
 * This example shows how vtkgdcmm can be connected to Kitware.VTK Activiz product.
 * Three (3) arguments are required:
 * 1. Input DICOM file (SWIG)
 * 2. Temporary PNG (intermediate) file (Activiz)
 * 3. Final DICOM file (SWIG)
 *
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/HelloActiviz.exe ~/Creatis/gdcmData/test.acr out.png toto.dcm
 *
 * Footnote:
 * this test originally used vtkBMPWriter / vtkBMPReader combination to store intermediate
 * image file, but BMP file are 24bits by default. Instead use PNG format which supports seems
 * to be closer to what was expected in this simple test.
 */
public class HelloActiviz
{
    // Does not work with Activiz.NET-5.4.0.455-Linux-x86_64-Personal
    /*
    static void ConnectSWIGToActiviz(Kitware.VTK.vtkImageExport imgin, Kitware.VTK.vtkImageImport imgout)
    {
        imgout.SetUpdateInformationCallback(imgin.GetUpdateInformationCallback());
        imgout.SetPipelineModifiedCallback(imgin.GetPipelineModifiedCallback());
        imgout.SetWholeExtentCallback(imgin.GetWholeExtentCallback());
        imgout.SetSpacingCallback(imgin.GetSpacingCallback());
        imgout.SetOriginCallback(imgin.GetOriginCallback());
        imgout.SetScalarTypeCallback(imgin.GetScalarTypeCallback());
        imgout.SetNumberOfComponentsCallback(imgin.GetNumberOfComponentsCallback());
        imgout.SetPropagateUpdateExtentCallback(imgin.GetPropagateUpdateExtentCallback());
        imgout.SetUpdateDataCallback(imgin.GetUpdateDataCallback());
        imgout.SetDataExtentCallback(imgin.GetDataExtentCallback());
        imgout.SetBufferPointerCallback(imgin.GetBufferPointerCallback());
        imgout.SetCallbackUserData(imgin.GetCallbackUserData());
    }
    */
}

```



```

static Kitware.VTK.vtkImageData ConnectSWIGToActiviz(vtkgdc.vtkImageData imgin)
{
    HandleRef rawCppThis = imgin.GetCppThis();
    Kitware.VTK.vtkImageData imgout = new Kitware.VTK.vtkImageData( rawCppThis.Handle, false, false);
    return imgout;
}

static vtkgdc.vtkImageData ConnectActivizToSWIG(Kitware.VTK.vtkImageData imgin)
{
    HandleRef rawCppThis = imgin.GetCppThis();
    vtkgdc.vtkImageData imgout = new vtkgdc.vtkImageData( rawCppThis );
    return imgout;
}

public static int Main(string[] args)
{
    string filename = args[0];
    string outfilename = args[1];

    // Step 1. Test SWIG -> Activiz
    vtkGDCMImageReader reader = vtkGDCMImageReader.
        New();
    reader.SetFileName( filename );
    //reader.Update(); // DO NOT call Update to check pipeline execution

    Kitware.VTK.vtkImageData imgout = ConnectSWIGToActiviz( reader.GetOutput() );

    System.Console.WriteLine( imgout.ToString() ); // not initialized as expected

    vtkPNGWriter writer = new vtkPNGWriter();
    writer.SetInput( imgout );
    writer.SetFileName( outfilename );
    writer.Write();

    // Step 2. Test Activiz -> SWIG
    vtkPNGReader bmpreader = new vtkPNGReader();
    bmpreader.SetFileName( outfilename );
    //bmpreader.Update(); // DO NOT update to check pipeline execution

    System.Console.WriteLine( bmpreader.GetOutput().ToString() ); // not initialized as expected

    vtkgdc.vtkImageData imgout2 = ConnectActivizToSWIG(bmpreader.GetOutput());

    System.Console.WriteLine( imgout2.ToString() ); // not initialized as expected

    Kitware.VTK.vtkMedicalImageProperties prop = new Kitware.VTK.vtkMedicalImageProperties();
    prop.SetModality( "MR" );

    string outfilename2 = args[2];
    vtkGDCMImageWriter writer2 = vtkGDCMImageWriter.
        New();
    writer2.SetMedicalImageProperties( prop.CastToActiviz() );
    writer2.SetFileName( outfilename2 );
    writer2.SetInput( imgout2 );
    writer2.Write();

    return 0;
}

```

## 29.86 HelloActiviz2.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdc.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====

```

```

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;

/*
 * Usage:
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/HelloActiviz2.exe gdcData/test.acr bla.png bla2.dcm
 */

/*
 * From the outside view, no-one can detect that object pass to/from
 * vtkGDCMImageWriter/vtkGDCMImageReader are not Activiz object.
 *
 * TODO: Test Command/Observer
 */
public class HelloActiviz2
{
    public static int Main(string[] args)
    {
        string filename = args[0];
        string outfilename = args[1];
        string outfilename2 = args[2];

        vtkGDCMImageReader reader = new Kitware.VTK.GDCM.
            vtkGDCMImageReader();
        reader.SetFileName( filename );

        // When calling multiple times creation of C# object from the same C++ object it triggers a:
        //error: potential refcounting error: Duplicate rawCppThis - weak reference that is still alive. Attempting
        //to add '0x00b2dc10' again.
        //    Allowing new wrapped object to take over table key...
        //    Original object should *not* have been destroyed while we still had it in our table without
        //    notifying us...
        //reader.GetOutput();
        //reader.GetOutput();

        System.Console.WriteLine( reader.ToString() ); // Test the ToString compat with Activiz

        vtkGDCMImageWriter writer = new vtkGDCMImageWriter();
        writer.SetInput( reader.GetOutput() );
        writer.SetFileName( outfilename2 );
        writer.Write();

        System.Console.WriteLine( reader.GetOutput().ToString() ); // Test the ToString compat with Activiz

        System.Console.WriteLine( writer.ToString() ); // Test the ToString compat with Activiz

        vtkPNGWriter pngwriter = new vtkPNGWriter();
        pngwriter.SetInput( reader.GetOutput() );
        pngwriter.SetFileName( outfilename );
        pngwriter.Write();

        // at that point the .Write() should have triggered an Update() on the reader:
        if( reader.GetImageFormat() == vtkgdc.VTK_LUMINANCE ) // MONOCHROME2
        {
            System.Console.WriteLine( "Image is MONOCHROME2" ); //
        }

        vtkPNGReader bmpreader = new vtkPNGReader();
        bmpreader.SetFileName( outfilename );

        vtkMedicalImageProperties prop = new vtkMedicalImageProperties();
        prop.SetModality( "MR" );

        vtkMatrix4x4 dircos = reader.GetDirectionCosines();
        dircos.Invert();

        vtkGDCMImageWriter writer2 = new vtkGDCMImageWriter();
        writer2.SetFileName( outfilename2 );
        writer2.SetDirectionCosines( dircos );
        writer2.SetMedicalImageProperties( prop );
        writer2.SetInput( bmpreader.GetOutput() );
        writer2.Write();

        return 0;
    }
}

```

## 29.87 HelloActiviz3.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;

/*
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/HelloActiviz3.exe ~/Creatis/gdcmData/test.acr
 */
public class HelloActiviz3
{
    public static int Main(string[] args)
    {
        string filename = args[0];

        vtkGDCMImageReader reader = vtkGDCMImageReader.
            New();
        vtkStringArray array = vtkStringArray.New();
        array.InsertNextValue(filename);

        reader.SetFileNames(array);
        reader.Update();

        //System.Console.WriteLine(reader.GetOutput());

        vtkRenderWindowInteractor iren = vtkRenderWindowInteractor.New();

        vtkImageViewer2 viewer = vtkImageViewer2.New();
        viewer.SetInput(reader.GetOutput());
        viewer.SetupInteractor(iren);
        viewer.SetSize(600, 600);
        viewer.Render();

        iren.Initialize();
        iren.Start();

        return 0;
    }
}

```

## 29.88 HelloActiviz4.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;

/*
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/HelloActiviz4.exe ~/Creatis/gdcmData/test.acr
 */
public class HelloActiviz4

```

```

{
    public static int Main(string[] args)
    {
        string filename = args[0];

        vtkGDCMImageReader reader = new vtkGDCMImageReader();
        vtkStringArray array = vtkStringArray.New();
        array.InsertNextValue(filename);

        reader.SetFileNames(array);
        reader.Update();

        //System.Console.Write(reader.GetOutput());

        vtkRenderWindowInteractor iren = vtkRenderWindowInteractor.New();

        vtkImageViewer viewer = vtkImageViewer.New();
        viewer.SetInput(reader.GetOutput());
        viewer.SetupInteractor(iren);
        viewer.SetSize(600, 600);
        viewer.Render();

        iren.Initialize();
        iren.Start();

        return 0;
    }
}

```

## 29.89 HelloActiviz5.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;

// The command line arguments are:
// -I      => run in interactive mode; unless this is used, the program will
//          not allow interaction and exit
// -D <path> => path to the data; the data should be in <path>/Data/

/*
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/HelloActiviz5.exe -I
 */
public class HelloActiviz5
{
    public static int Main(string[] args)
    {
        vtkTesting testHelper = vtkTesting.New();
        for (int cc = 0; cc < args.Length; cc++)
        {
            //testHelper.AddArguments(argc, const_cast<const char **>(argv));
            //System.Console.Write( "args: " + args[cc] + "\n" );
            testHelper.AddArgument( args[cc] );
        }
        if ( testHelper.IsFlagSpecified("-D") != 0 )
        {
            string VTK_DATA_ROOT = vtkGDCMTesting.GetVTKDataRoot();
            if ( VTK_DATA_ROOT != null )
            {
                //System.Console.Write( "VTK_DATA_ROOT: " + VTK_DATA_ROOT + "\n" );
                testHelper.SetDataRoot(VTK_DATA_ROOT);
                testHelper.AddArgument("-D");
                testHelper.AddArgument(VTK_DATA_ROOT);
            }
        }
    }
}

```

```

    }

    string dataRoot = testHelper.GetDataRoot();
    string filename = dataRoot;
    filename += "/Data/mr.001";

    vtkDirectory dir = vtkDirectory.New();
    if( dir.FileIsDirectory( dataRoot ) == 0 )
    {
        filename = vtkGDCMTesting.GetGDCMDataRoot() + "/test.acr";
    }
    //System.Console.Write( "dataRoot: " + dataRoot + "\n" );
    System.Console.Write( "filename being used is: " + filename + "\n" );

    vtkGDCMImageReader reader = vtkGDCMImageReader.
        New();
    vtkStringArray array = vtkStringArray.New();
    array.InsertNextValue(filename);
    reader.SetFileNames(array);
    reader.Update();

    System.Console.Write(reader.GetOutput());

    vtkRenderWindowInteractor iren = vtkRenderWindowInteractor.New();

    vtkRenderer ren1 = vtkRenderer.New();
    vtkRenderWindow renWin = vtkRenderWindow.New();
    renWin.AddRenderer(ren1);

    vtkImageActor actor = vtkImageActor.New();

    vtkImageMapToWindowLevelColors coronalColors = vtkImageMapToWindowLevelColors.
        New();
    coronalColors.SetInput(reader.GetOutput());

    actor.SetInput(coronalColors.GetOutput());

    ren1.AddActor(actor);
    iren.SetRenderWindow(renWin);

    iren.Initialize();

    renWin.Render();

    int retVal = testHelper.IsInteractiveModeSpecified();

    if( retVal != 0 )
    {
        iren.Start();
    }

    return 0;
}

```

## 29.90 HelloSimple.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

/*
 * Compilation:
 * $ CLASSPATH=gdcm.jar javac ../../gdcm/Examples/Java/HelloSimple.java -d .
 *
 * Usage:
 * $ LD_LIBRARY_PATH=. CLASSPATH=gdcm.jar:. java HelloSimple gdcmData/012345.002.050.dcm
 */

```

```

import gdc.*;

public class HelloSimple
{
    public static void main(String[] args) throws Exception
    {
        String filename = args[0];
        Reader reader = new Reader();
        reader.SetFileName( filename );
        boolean ret = reader.Read();
        if( !ret )
        {
            throw new Exception("Could not read: " + filename );
        }
        File f = reader.GetFile();
        DataSet ds = f.GetDataSet();

        System.out.println( ds.toString() );

        System.out.println("Success reading: " + filename );
    }
}

```

## 29.91 HelloVizWorld.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdc.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

/*
 * Basic example for dealing with a DICOM file that contains an Image
 * (read: Pixel Data element)
 */

#include "gdcImageReader.h"
#include "gdcImageWriter.h"
#include "gdcImage.h"
#include "gdcPhotometricInterpretation.h"

#include <iostream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    // Instantiate the image reader:
    gdc::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }
    // If we reach here, we know for sure 2 things:
    // 1. It is a valid DICOM
    // 2. And it contains an Image !

    // The output of superclass gdc::Reader is a gdc::File
    //gdc::File &file = reader.GetFile();

    // The other output of gdc::ImageReader is a gdc::Image
    const gdc::Image &image = reader.GetImage();

```

```

// Let's get some property from the image:
unsigned int ndim = image.GetNumberOfDimensions();
// Dimensions of the image:
const unsigned int *dims = image.GetDimensions();
// Origin
const double *origin = image.GetOrigin();
const gdcm::PhotometricInterpretation &pi = image.
    GetPhotometricInterpretation();
for(unsigned int i = 0; i < ndim; ++i)
{
    std::cout << "Dim(" << i << "): " << dims[i] << std::endl;
}
for(unsigned int i = 0; i < ndim; ++i)
{
    std::cout << "Origin(" << i << "): " << origin[i] << std::endl;
}
std::cout << "PhotometricInterpretation: " << pi << std::endl;

// Write the modified DataSet back to disk
gdcm::ImageWriter writer;
writer.SetImage( image );
writer.SetFileName( outfilename );
//writer.SetFile( file ); // We purposely NOT copy the meta information from the input
// file, and instead only pass the image
if( !writer.Write() )
{
    std::cerr << "Could not write: " << outfilename << std::endl;
    return 1;
}

return 0;
}

```

## 29.92 HelloVTKWorld.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
using vtkgdcm;

/*
 * This test only test the SWIG/VTK part, you do not need Activiz
 */
public class HelloVTKWorld
{
    public static int Main(string[] args)
    {
        {
            string filename = args[0];
            vtkGDCMImageReader reader = vtkGDCMImageReader.
                New();
            reader.SetFileName( filename );
            reader.Update();

            vtkMedicalImageProperties prop = reader.GetMedicalImageProperties();
            System.Console.WriteLine( prop.GetPatientName() ); //

            if( reader.GetImageFormat() == vtkgdcm.vtkgdcm.VTK_LUMINANCE ) // MONOCHROME2
            {
                System.Console.WriteLine( "Image is MONOCHROME2" ); //
            }

            // Just for fun, invert the direction cosines, output should reflect that:
            vtkMatrix4x4 dircos = reader.GetDirectionCosines();
            dircos.Invert();

            string outfilename = args[1];

```

```

    vtkGDCMImageWriter writer = vtkGDCMImageWriter.
        New();
    writer.SetMedicalImageProperties( reader.GetMedicalImageProperties() );
    writer.SetDirectionCosines( dircos );
    writer.SetShift( reader.GetShift() );
    writer.SetScale( reader.GetScale() );
    writer.SetImageFormat( reader.GetImageFormat() );
    writer.SetFileName( outfilename );
    //writer.SetInputConnection( reader.GetOutputPort() ); // new
    writer.SetInput( reader.GetOutput() ); // old
    writer.Write();

    return 0;
}

```

## 29.93 HelloVTKWorld.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
// We are required to call the package 'vtk' eventhough I (MM) would have preferred
// an import statement along the line of:
// import vtkgdc.*;
import vtk.*;

/*
 * Compilation:
 * CLASSPATH=vtkgdc.jar:/usr/share/java/vtk.jar javac HelloVTKWorld.java
 *
 * Usage:
 * LD_LIBRARY_PATH=/usr/lib/jvm/java-6-openjdk/jre/lib/amd64/xawt:/usr/lib/jni:. CLASSPATH=/usr/share/java/
 *   vtk.jar:vtkgdc.jar:gdcm.jar:. java HelloVTKWorld gdcmData/012345.002.050.dcm bla.dcm
 */
public class HelloVTKWorld
{
    static {
        System.loadLibrary("vtkCommonJava");
        System.loadLibrary("vtkFilteringJava");
        System.loadLibrary("vtkIOJava");
        System.loadLibrary("vtkImagingJava");
        System.loadLibrary("vtkGraphicsJava");
        System.loadLibrary("vtkgdcJava");
        try {
            System.loadLibrary("vtkRenderingJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkHybrid, skipping...");
        }
        try {
            System.loadLibrary("vtkHybridJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkHybrid, skipping...");
        }
        try {
            System.loadLibrary("vtkVolumeRenderingJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkVolumeRendering, skipping...");
        }
    }

    public static void main(String[] args)
    {
        String filename = args[0];
        vtkGDCMImageReader reader = new vtkGDCMImageReader();
        reader.SetFileName( filename );
        reader.Update();
    }
}

```



```

        vtkMedicalImageProperties prop = reader.GetMedicalImageProperties();
        System.out.println( prop.GetPatientName() ); //

//      if( reader.GetImageFormat() == vtkgdc.VtkGDCM.VTK_LUMINANCE ) // MONOCHROME2
//      {
//          System.out.println( "Image is MONOCHROME2" ); //
//      }

// Just for fun, invert the direction cosines, output should reflect that:
vtkMatrix4x4 dircos = reader.GetDirectionCosines();
dircos.Invert();

// We need to maintain in sync information stored in vtkMedicalImageProperties:
double[] cosines = new double[6];
cosines[0] = dircos.GetElement(0,0);
cosines[1] = dircos.GetElement(1,0);
cosines[2] = dircos.GetElement(2,0);
cosines[3] = dircos.GetElement(0,1);
cosines[4] = dircos.GetElement(1,1);
cosines[5] = dircos.GetElement(2,1);
reader.GetMedicalImageProperties().SetDirectionCosine( cosines );

String outfilename = args[1];
vtkGDCMImageWriter writer = new vtkGDCMImageWriter();
writer.SetMedicalImageProperties( reader.GetMedicalImageProperties() );
writer.SetDirectionCosines( dircos );
writer.SetShift( reader.GetShift() );
writer.SetScale( reader.GetScale() );
writer.SetImageFormat( reader.GetImageFormat() );
writer.SetFileName( outfilename );
//writer.SetInputConnection( reader.GetOutputPort() ); // new
writer.SetInput( reader.GetOutput() ); // old
writer.Write();

System.out.println("Success reading: " + filename );
}
}

```

## 29.94 HelloVTKWorld2.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using vtkgdc;

/*
 * This test only test the SWIG/VTK part, you do not need Activiz
 */
public class HelloVTKWorld2
{
    public static int Main(string[] args)
    {
        string VTK_DATA_ROOT = vtkGDCMTesting.GetVTKDataRoot();

        vtkVoxel6Reader reader = vtkVoxel6Reader.New();
        reader.SetDataDimensions(64, 64);
        reader.SetDataByteOrderToLittleEndian();
        reader.SetFilePrefix(VTK_DATA_ROOT + "/Data/headsq/quarter");
        reader.SetImageRange(1, 93);
        reader.SetDataSpacing(3.2, 3.2, 1.5);

        vtkImageCast cast = vtkImageCast.New();
        cast.SetInput( reader.GetOutput() );
        cast.SetOutputScalarTypeToUnsignedChar();

        // By default this is creating a Multiframe Grayscale Word Secondary Capture Image Storage
    }
}

```

```

    vtkGDCMImageWriter writer = vtkGDCMImageWriter.
        New();
    writer.SetFileName( "headsq.dcm" );
    writer.SetInput( reader.GetOutput() );
    // cast -> Multiframe Grayscale Byte Secondary Capture Image Storage
    // writer.SetInput( cast.GetOutput() );
    writer.SetFileDimensionality( 3 );
    writer.Write();

    return 0;
}

```

## 29.95 HelloWorld.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
/*
 * This example is ... guess what this is for :)
 */

#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmAttribute.h"

#include <iostream>

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];

    // Instantiate the reader:
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }

    // If we reach here, we know for sure only 1 thing:
    // It is a valid DICOM file (potentially an old ACR-NEMA 1.0/2.0 file)
    // (Maybe, it's NOT a Dicom image -could be a DICOMDIR, a RTSTRUCT, etc-)

    // The output of gdcm::Reader is a gdcm::File
    gdcm::File &file = reader.GetFile();

    // the dataset is the the set of element we are interested in:
    gdcm::DataSet &ds = file.GetDataSet();

    // Construct a static(*) type for Image Comments :
    gdcm::Attribute<0x0020,0x4000> imagecomments;
    imagecomments.SetValue( "Hello, World !" );

    // Now replace the Image Comments from the dataset with our:
    ds.Replace( imagecomments.GetAsDataElement() );

    // Write the modified DataSet back to disk
    gdcm::Writer writer;
    writer.CheckFileMetaInformationOff(); // Do not attempt to reconstruct the

```

```

        file meta to preserve the file
        // as close to the original as possible.
writer.SetFileName( outfilename );
writer.SetFile( file );
if( !writer.Write() )
{
    std::cerr << "Could not write: " << outfilename << std::endl;
    return 1;
}

return 0;
}

/*
 * (*) static type, means that extra DICOM information VR & VM are computed at compilation time.
 * The compiler is deducing those values from the template arguments of the class.
 */

```

## 29.96 HelloWorld.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Hello World !
17 """
18
19 import gdcm
20 import sys
21
22 if __name__ == "__main__":
23
24     # verbosity:
25     #gdcm.Trace.DebugOn()
26     #gdcm.Trace.WarningOn()
27     #gdcm.Trace.ErrorOn()
28
29     # Get the filename from the command line
30     filename = sys.argv[1]
31
32     # Instanciate a gdcm.Reader
33     # This is the main class to handle any type of DICOM object
34     # You should check for gdcm.ImageReader for reading specifically DICOM Image file
35     r = gdcm.Reader()
36     r.SetFileName( filename )
37     # If the reader fails to read the file, we should stop !
38     if not r.Read():
39         print "Not a valid DICOM file"
40         sys.exit(1)
41
42     # Get the DICOM File structure
43     file = r.GetFile()
44
45     # Get the DataSet part of the file
46     dataset = file.GetDataSet()
47
48     # Ok let's print it !
49     print dataset
50
51     # Use StringFilter to print a particular Tag:
52     sf = gdcm.StringFilter()
53     sf.SetFile(r.GetFile())
54
55     # Check if Attribute exist
56     print dataset.FindElement( gdcm.Tag(0x0028,0x0010) )
57

```

```

58  # Let's print it as string pair:
59  print sf.ToStringPair(gdcm.Tag(0x0028,0x0010))

```

## 29.97 iU22tomultisc.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * iU22 Raw Data extractor
 */
#include "gdcmReader.h"
#include "gdcmImageWriter.h"
#include "gdcmAttribute.h"
#include "gdcmPrivateTag.h"

#include <math.h>

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    // IM_001
    const char *filename = argv[1];

    gdcm::Reader reader; // Do not use ImageReader
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }

    // * The data is simply 8-bit unsigned in the obvious x/y/z order
    // * 200D,300B contains the data
    // * 200D,3001 contains the no. of voxels (416,412,256 in this case)
    // * 200D,3003 contains the voxel sizes (0.156184527398215 /
    // 0.1223749613981957 / 0.328479990704639 in this case)

    const gdcm::File &file = reader.GetFile();
    const gdcm::DataSet &ds = file.GetDataSet();
    const gdcm::PrivateTag trawdataus( 0x200d, 0x0b, "Philips US Imaging DD 033" );
    const gdcm::DataElement &rawdataus = ds.GetDataElement( trawdataus );

    const gdcm::PrivateTag tcolsrowsframes( 0x200d, 0x01, "Philips US Imaging DD 036" );
    const gdcm::DataElement &colsrowsframes = ds.GetDataElement(
        tcolsrowsframes );
    // const gdcm::PrivateTag tcolsrowsframes( 0x200d, 0x02, "Philips US Imaging DD 036" );
    // this is just a duplicate previous tag.
    const gdcm::PrivateTag tvoxelspacing( 0x200d, 0x03, "Philips US Imaging DD 036" );
    const gdcm::DataElement &voxelspacing = ds.GetDataElement( tvoxelspacing )
    ;

    gdcm::Element<gdcm::VR::DS, gdcm::VM::VM3> dims; // Use DS to
        interpret value stored in LO
    dims.SetFromDataElement( colsrowsframes );

    gdcm::Element<gdcm::VR::DS, gdcm::VM::VM3> spacing;
    spacing.SetFromDataElement( voxelspacing );

    gdcm::ImageWriter writer;

    gdcm::Image &image = writer.GetImage();
    image.SetNumberOfDimensions( 3 ); // good default
    image.SetDimension(0, (unsigned int)dims[0] );
    image.SetDimension(1, (unsigned int)dims[1] );
    image.SetDimension(2, (unsigned int)dims[2] );
    image.SetSpacing(0, spacing[0] );

```

```

image.SetSpacing(1, spacing[1] );
image.SetSpacing(2, spacing[2] );
gdcmm::PixelFormat pixeltype = gdcmm::PixelFormat::UINT8;

gdcmm::PhotometricInterpretation pi;
pi = gdcmm::PhotometricInterpretation::MONOCHROME2;
image.SetPhotometricInterpretation( pi );
image.SetPixelFormat( pixeltype );

image.SetDataElement( rawdataus );

std::string outfilename = "outiu22.dcm";

gdcmm::DataElement de( gdcmm::Tag(0x8,0x16) ); // SOP Class UID
de.SetVR( gdcmm::VR::UI );
gdcmm::MediaStorage ms(
    gdcmm::MediaStorage::UltrasoundMultiFrameImageStorage
);
// gdcmm::MediaStorage::MultiframeGrayscaleByteSecondaryCaptureImageStorage );
de.SetByteValue( ms.GetString(), (uint32_t)strlen(ms.
    GetString()));
writer.GetFile().GetDataSet().Replace( de );

writer.SetFileName( outfilename.c_str() );
if( !writer.Write() )
{
    std::cerr << "could not write: " << outfilename << std::endl;
    return 1;
}

return 0;
}

```

## 29.98 LargeVRDSExplicit.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmmReader.h"
#include "gdcmmWriter.h"
#include "gdcmmAttribute.h"
#include "gdcmmFileExplicitFilter.h"
#include "gdcmmSequenceOfItems.h"

bool interpolate(const double * pts, size_t npts, std::vector<double> &out )
{
    out.clear();
    for(size_t i = 0; i < 2*npts; ++i )
    {
        const size_t j = i / 2;
        if( i % 2 )
        {
            if( j != npts - 1 )
            {
                assert( 3*j+5 < 3*npts );
                const double midpointx = (pts[3*j+0] + pts[3*j+3]) / 2;
                const double midpointy = (pts[3*j+1] + pts[3*j+4]) / 2;
                const double midpointz = (pts[3*j+2] + pts[3*j+5]) / 2;
                out.push_back( midpointx );
                out.push_back( midpointy );
                out.push_back( midpointz );
            }
        }
        else
        {
            assert( j < npts );

```

```

        out.push_back( pts[3*j+0] );
        out.push_back( pts[3*j+1] );
        out.push_back( pts[3*j+2] );
    }
}
assert( out.size() == 2 * npts * 3 - 3 );
return true;
}

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    const char *outfilename = argv[2];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();

    gdcm::FileExplicitFilter fef;
    //fef.SetChangePrivateTags( changeprivatetags );
    fef.SetFile( reader.GetFile() );
    if( !fef.Change() )
    {
        std::cerr << "Failed to change: " << filename << std::endl;
        return 1;
    }

    // (3006,0039) SQ (Sequence with undefined length #=4)      # u/l, 1 ROIContourSequence
    gdcm::Tag tag(0x3006,0x0039);

    const gdcm::DataElement &roicsq = ds.GetDataElement( tag );
    gdcm::SmartPointer<gdcm::SequenceOfItems> sqi = roicsq.
        GetValueAsSQ();
    //sqi->SetNumberOfItems( 1 );
    const gdcm::Item & item = sqi->GetItem(1); // Item start at #1
    const gdcm::DataSet& nestedds = item.GetNestedDataSet();

    gdcm::Tag tcsq(0x3006,0x0040);
    if( !nestedds.FindDataElement( tcsq ) )
    {
        return 0;
    }
    const gdcm::DataElement& csq = nestedds.GetDataElement( tcsq );
    gdcm::SmartPointer<gdcm::SequenceOfItems> sqi2 = csq.
        GetValueAsSQ();
    if( !sqi2 || !sqi2->GetNumberOfItems() )
    {
        return 0;
    }
    //unsigned int nitems = sqi2->GetNumberOfItems();
    gdcm::Item & item2 = sqi2->GetItem(1); // Item start at #1

    gdcm::DataSet& nestedds2 = item2.GetNestedDataSet();
    //item2.SetVLToUndefined();
    //std::cout << nestedds2 << std::endl;
    // (3006,0050) DS [43.57636\65.52504\10.0\46.043102\62.564945\10.0\49.126537\60.714... # 398,48
    ContourData
    gdcm::Tag tcontourdata(0x3006,0x0050);
    const gdcm::DataElement & contourdata = nestedds2.
        GetDataElement( tcontourdata );
    //std::cout << contourdata << std::endl;

    //const gdcm::ByteValue *bv = contourdata.GetByteValue();
    gdcm::Attribute<0x3006,0x0046> ncontourpoints;
    ncontourpoints.Set( nestedds2 );

    gdcm::Attribute<0x3006,0x0050> at;
    at.SetFromDataElement( contourdata );
    const double* pts = at.GetValues();
    unsigned int npts = at.GetNumberOfValues() / 3;

```

```

std::vector<double> out( pts, pts + npts * 3 );
std::vector<double> out2;

//const unsigned int niter = 7;
const unsigned int niter = 8;
for( unsigned int i = 0; i < niter; ++i)
{
    //bool b =
    interpolate(&out[0], out.size() / 3, out2);
    //const double *pout = &out[0];
    out = out2;
    out2.clear();
}
assert( out.size() % 3 == 0 );

gdcm::Attribute<0x3006,0x0050> at_interpolate;
at_interpolate.SetNumberOfValues( (unsigned int)(out.size() / 3) );
at_interpolate.SetValues( &out[0], (uint32_t)out.size() );

ncontourpoints.SetValue( at_interpolate.GetNumberOfValues() / 3 );
nestedds2.Replace( at_interpolate.GetAsDataElement() );
nestedds2.Replace( ncontourpoints.GetAsDataElement() );

//assert(0);

// Let's take item one and subdivide it

gdcm::TransferSyntax ts =
    gdcm::TransferSyntax::ImplicitVRLittleEndian;
ts = gdcm::TransferSyntax::ExplicitVRLittleEndian;

gdcm::FileMetaInformation &fmi = file.GetHeader();
const char *tsuid = gdcm::TransferSyntax::GetTSString( ts );
// const char * is ok since padding is \0 anyway...
gdcm::DataElement de( gdcm::Tag(0x0002,0x0010) );
de.SetByteValue( tsuid, (uint32_t)strlen(tsuid) );
de.SetVR( gdcm::Attribute<0x0002, 0x0010>::GetVR() );
fmi.Replace( de );
fmi.Remove( gdcm::Tag(0x0002,0x0012) ); // will be regenerated
fmi.Remove( gdcm::Tag(0x0002,0x0013) ); // ' ' ' '
fmi.SetDataSetTransferSyntax(ts);

gdcm::Writer w;
w.SetFile( file );
w.SetFileName( outfilename );
if (!w.Write() )
{
    return 1;
}

return 0;
}

```

## 29.99 MagnifyFile.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkGDCMImageWriter.h"
#include "vtkImageData.h"
#include "vtkImageMagnify.h"
#include "vtkImageCast.h"

#include "gdcmTesting.h"
#include "gdcmSystem.h"

```

```
// This is a simple test to magnify an image that is known to give excellent
// compression ratio. This will be our test for those large image
int main(int, char *[])
{
    const char *directory = gdcm::Testing::GetDataRoot();
    if(!directory) return 1;
    std::string file = std::string(directory) + "/test.acr";
    std::cout << file << std::endl;
    if( !gdcm::System::FileExists( file.c_str() ) ) return 1;

    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( file.c_str() );
    reader->Update();
    //reader->GetOutput()->Print( std::cout );

    vtkImageCast *cast = vtkImageCast::New();
    #if (VTK_MAJOR_VERSION >= 6)
        cast->SetInputConnection( reader->GetOutputPort() );
    #else
        cast->SetInput( reader->GetOutput() );
    #endif
    cast->SetOutputScalarTypeToUnsignedShort();

    vtkImageMagnify *magnify = vtkImageMagnify::New();
    #if (VTK_MAJOR_VERSION >= 6)
        magnify->SetInputConnection( cast->GetOutputPort() );
    #else
        magnify->SetInput( cast->GetOutput() );
    #endif
    magnify->SetInterpolate( 1 );
    magnify->SetInterpolate( 0 );
    int factor = 100;
    magnify->SetMagnificationFactors (factor, factor, 1);

    vtkGDCMImageWriter *writer = vtkGDCMImageWriter::New();
    writer->SetFileName( "/tmp/bla.dcm" );
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputConnection( magnify->GetOutputPort() );
    #else
        writer->SetInput( magnify->GetOutput() );
    #endif
    writer->SetImageFormat( reader->GetImageFormat() );
    writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
    writer->SetDirectionCosines( reader->GetDirectionCosines() );
    writer->SetShift( reader->GetShift() );
    writer->SetScale( reader->GetScale() );
    writer->Write();

    // TODO:
    //vtkImageAppendComponents.h

    reader->Delete();
    magnify->Delete();
    writer->Delete();

    return 0;
}
```

## 29.100 ManipulateFile.cs

```
/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
```



```

* $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
* $ mono bin/ManipulateFile.exe gdcmData/012345.002.050.dcm out.dcm
*/
using System;
using gdcm;

public class ManipulateFile
{
    public static int Main(string[] args)
    {
        string file1 = args[0];
        string file2 = args[1];
        Reader reader = new Reader();
        reader.SetFileName( file1 );
        bool ret = reader.Read();
        if( !ret )
        {
            return 1;
        }

        Anonymizer ano = new Anonymizer();
        ano.SetFile( reader.GetFile() );
        ano.RemovePrivateTags();
        ano.RemoveGroupLength();
        Tag t = new Tag(0x10,0x10);
        ano.Replace( t, "GDCM^Csharp^Test^Hello^World" );

        UIDGenerator g = new UIDGenerator();
        ano.Replace( new Tag(0x0008,0x0018), g.Generate() );
        ano.Replace( new Tag(0x0020,0x000d), g.Generate() );
        ano.Replace( new Tag(0x0020,0x000e), g.Generate() );
        ano.Replace( new Tag(0x0020,0x0052), g.Generate() );

        Writer writer = new Writer();
        writer.SetFileName( file2 );
        writer.SetFile( ano.GetFile() );
        ret = writer.Write();
        if( !ret )
        {
            return 1;
        }

        return 0;
    }
}

```

## 29.101 ManipulateFile.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python ManipulateFile.py input.dcm output.dcm
19
20 Footnote:
21 GDCM 1.2.x would create incorrect Multiframe MR Image Storage file. Try to recover from
22 the issues to recreate a MultiframeGrayscaleByteSecondaryCaptureImageStorage file.
23 e.g:
24
25 python ManipulateFile.py Insight/Testing/Temporary/itkGDCMImageIOTest5-j2k.dcm manipulated.dcm
26 """
27
28 import sys
29 import gdcm

```

```

30
31 if __name__ == "__main__":
32
33     file1 = sys.argv[1]
34     file2 = sys.argv[2]
35
36     r = gdcm.Reader()
37     r.SetFileName( file1 )
38     if not r.Read():
39         sys.exit(1)
40
41     ano = gdcm.Anonymizer()
42     ano.SetFile( r.GetFile() )
43     ano.RemovePrivateTags()
44     ano.Remove( gdcm.Tag(0x0032,0x1030) )
45     ano.Remove( gdcm.Tag(0x008,0x14) )
46     ano.Remove( gdcm.Tag(0x008,0x1111) )
47     ano.Remove( gdcm.Tag(0x008,0x1120) )
48     ano.Remove( gdcm.Tag(0x008,0x1140) )
49     ano.Remove( gdcm.Tag(0x10,0x21b0) )
50     ano.Empty( gdcm.Tag(0x10,0x10) )
51     ano.Empty( gdcm.Tag(0x10,0x20) )
52     ano.Empty( gdcm.Tag(0x10,0x30) )
53     ano.Empty( gdcm.Tag(0x20,0x10) )
54     ano.Empty( gdcm.Tag(0x32,0x1032) )
55     ano.Empty( gdcm.Tag(0x32,0x1033) )
56     ano.Empty( gdcm.Tag(0x40,0x241) )
57     ano.Empty( gdcm.Tag(0x40,0x254) )
58     ano.Empty( gdcm.Tag(0x40,0x253) )
59     ano.Empty( gdcm.Tag(0x40,0x1001) )
60     ano.Empty( gdcm.Tag(0x8,0x80) )
61     ano.Empty( gdcm.Tag(0x8,0x50) )
62     ano.Empty( gdcm.Tag(0x8,0x1030) )
63     ano.Empty( gdcm.Tag(0x8,0x103e) )
64     ano.Empty( gdcm.Tag(0x18,0x1030) )
65     ano.Empty( gdcm.Tag(0x38,0x300) )
66     g = gdcm.UIDGenerator()
67     ano.Replace( gdcm.Tag(0x0008,0x0018), g.Generate() )
68     ano.Replace( gdcm.Tag(0x0020,0x00d), g.Generate() )
69     ano.Replace( gdcm.Tag(0x0020,0x00e), g.Generate() )
70     ano.Replace( gdcm.Tag(0x0020,0x052), g.Generate() )
71     #ano.Replace( gdcm.Tag(0x0008,0x0016), "1.2.840.10008.5.1.4.1.1.7.2" )
72     ""
73     ano.Remove( gdcm.Tag(0x0018,0x0020) ) # ScanningSequence
74     ano.Remove( gdcm.Tag(0x0018,0x0021) ) # SequenceVariant
75     ano.Remove( gdcm.Tag(0x0018,0x0022) ) # ScanOptions
76     ano.Remove( gdcm.Tag(0x0018,0x0023) ) # MRAcquisitionType
77     ano.Remove( gdcm.Tag(0x0018,0x0050) ) # SliceThickness
78     ano.Remove( gdcm.Tag(0x0018,0x0080) ) # RepetitionTime
79     ano.Remove( gdcm.Tag(0x0018,0x0081) ) # EchoTime
80     ano.Remove( gdcm.Tag(0x0018,0x0088) ) # SpacingBetweenSlices
81     ano.Remove( gdcm.Tag(0x0018,0x0091) ) # EchoTrainLength
82     ano.Remove( gdcm.Tag(0x0018,0x1164) ) # ImagerPixelSpacing
83
84     ano.Remove( gdcm.Tag(0x0020,0x0032) ) # Image Position (Patient)
85     ano.Remove( gdcm.Tag(0x0020,0x0037) ) # Image Orientation (Patient)
86     ano.Remove( gdcm.Tag(0x0020,0x0052) ) # Frame of Reference UID
87     ano.Remove( gdcm.Tag(0x0020,0x1040) ) # Position Reference Indicator
88
89     ano.Replace( gdcm.Tag(0x0028,0x0301), "NO" ) # Burned In Annotation
90
91     ano.Empty( gdcm.Tag(0x0020,0x0020) )
92
93     ano.Remove( gdcm.Tag(0x7fe0,0x0000) )
94
95     #ano.Empty( gdcm.Tag(0x0028,0x0009) ) # Frame Increment Pointer
96
97     #ano.Empty( gdcm.Tag(0x0028,0x1052) ) #<entry group="0028" element="1052" vr="DS" vm="1" name="Rescale
98     Intercept"/>
99     #ano.Empty( gdcm.Tag(0x0028,0x1053) ) #<entry group="0028" element="1053" vr="DS" vm="1" name="Rescale
100     Slope"/>
101     #ano.Replace( gdcm.Tag(0x0028,0x1054), "US" ) #<entry group="0028" element="1054" vr="LO" vm="1" name="
102     Rescale Type"/>
103
104     ano.Replace( gdcm.Tag(0x2050, 0x0020), "IDENTITY")
105     ""
106
107     w = gdcm.Writer()
108     w.SetFile( ano.GetFile() )
109     w.SetFileName( file2 )
110     if not w.Write():

```

```
108     sys.exit(1)
```

## 29.102 ManipulateSequence.py

```
1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python ManipulateSequence.py input.dcm output.dcm
19
20 This was tested using:
21
22 python ManipulateSequence.py gdcmData/D_CLUNIE_CT1_J2KI.dcm myoutput.dcm
23
24 This is a dummy example on how to modify a value set in a nested-nested dataset
25
26 WARNING:
27 Do not use as-is in production, this is just an example
28 This example works in an undefined length Item only (you need to explicitly recompute the length
   otherwise)
29 """
30
31 import sys
32 import gdcm
33
34 if __name__ == "__main__":
35
36     file1 = sys.argv[1]
37     file2 = sys.argv[2]
38
39     r = gdcm.Reader()
40     r.SetFileName( file1 )
41     if not r.Read():
42         sys.exit(1)
43
44     f = r.GetFile()
45     ds = f.GetDataSet()
46     tsis = gdcm.Tag(0x0008,0x2112) # SourceImageSequence
47     if ds.FindDataElement( tsis ):
48         sis = ds.GetDataElement( tsis )
49         #sqsis = sis.GetSequenceOfItems()
50         # GetValueAsSQ handle more cases
51         sqsis = sis.GetValueAsSQ()
52         if sqsis.GetNumberOfItems():
53             item1 = sqsis.GetItem(1)
54             nestedds = item1.GetNestedDataSet()
55             tprcs = gdcm.Tag(0x0040,0xa170) # PurposeOfReferenceCodeSequence
56             if nestedds.FindDataElement( tprcs ):
57                 prcs = nestedds.GetDataElement( tprcs )
58                 sqprcs = prcs.GetSequenceOfItems()
59                 if sqprcs.GetNumberOfItems():
60                     item2 = sqprcs.GetItem(1)
61                     nestedds2 = item2.GetNestedDataSet()
62                     # (0008,0104) LO [Uncompressed predecessor] # 24, 1 CodeMeaning
63                     tcm = gdcm.Tag(0x0008,0x0104)
64                     if nestedds2.FindDataElement( tcm ):
65                         cm = nestedds2.GetDataElement( tcm )
66                         mystr = "GDCM was here"
67                         cm.SetByteValue( mystr, gdcm.VL( len(mystr) ) )
68
69     w = gdcm.Writer()
70     w.SetFile( f )
71     w.SetFileName( file2 )
72     if not w.Write():
```

```
73     sys.exit(1)
```

## 29.103 MergeFile.py

```
1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #       This software is distributed WITHOUT ANY WARRANTY; without even
10 #       the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #       PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python MergeFile.py input1.dcm input2.dcm
19
20 It will produce a 'merge.dcm' output file, which contains all meta information from input1.dcm
21 and copy the Stored Pixel values from input2.dcm
22 This script even works when input2.dcm is a Secondary Capture and does not contains information
23 such as IOP and IPP...
24 """
25
26 import sys
27 import gdcm
28
29 if __name__ == "__main__":
30
31     file1 = sys.argv[1]
32     file2 = sys.argv[2]
33
34     r1 = gdcm.ImageReader()
35     r1.SetFileName( file1 )
36     if not r1.Read():
37         sys.exit(1)
38
39     r2 = gdcm.ImageReader()
40     r2.SetFileName( file2 )
41     if not r2.Read():
42         sys.exit(1)
43
44     # Image from r2 could be Secondary Capture and thus would not contains neither IPP nor IOP
45     # Instead always prefer to only copy the Raw Data Element.
46     # Warning ! Image need to be identical ! Only the value of Stored Pixel can be different.
47     r1.GetImage().SetDataElement( r2.GetImage().GetDataElement() )
48
49     w = gdcm.ImageWriter()
50     w.SetFile( r1.GetFile() )
51     #w.SetImage( r2.GetImage() ) # See comment above
52     w.SetImage( r1.GetImage() )
53
54     w.SetFileName( "merge.dcm" )
55     if not w.Write():
56         sys.exit(1)
57
58     sys.exit(0)
```

## 29.104 MergeTwoFiles.cxx

```
/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
```

```

    This software is distributed WITHOUT ANY WARRANTY; without even
    the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
    PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example will show how one can read in two DICOM files, use the dataset
 * from file1 and use image from file2 to save it in a 3rd file.
 *
 * Eg:
 * MergeTwoFiles gdcmlData/012345.002.050.dcm gdcmlData/test.acr merge.dcm
 */

#include "gdcmReader.h"
#include "gdcmImageReader.h"
#include "gdcmImageWriter.h"
#include "gdcmWriter.h"
#include "gdcmDataSet.h"
#include "gdcmAttribute.h"

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        return 1;
    }
    const char *file1 = argv[1];
    const char *file2 = argv[2];
    const char *file3 = argv[3];

    // Read file1
    gdcm::ImageReader reader1;
    reader1.SetFileName( file1 );
    if( !reader1.Read() )
    {
        return 1;
    }

    // Read file2
    gdcm::ImageReader reader2;
    reader2.SetFileName( file2 );
    if( !reader2.Read() )
    {
        return 1;
    }

    // Ok now let's take the DataSet from file1 and the Image from file2
    // Warning: if file2 is -for example- a Secondary Capture Storage, then it has no
    // Image Orientation (Patient) thus any Image Orientation (Patient) from file1
    // will be discarded...

    // let's be fancy. In case reader2 contains explicit, but reader1 is implicit
    // we would rather see an implicit output
    if( reader1.GetFile().GetHeader().GetDataSetTransferSyntax() ==
        gdcm::TransferSyntax::ImplicitVRLittleEndian )
    {
        reader2.GetImage().SetTransferSyntax(
            gdcm::TransferSyntax::ImplicitVRLittleEndian );
    }

    gdcm::ImageWriter writer;
    writer.SetFileName( file3 );
    writer.SetFile( reader1.GetFile() );
    // ImageWriter will always use all of gdcm::Image information and override anything wrong from
    // reader1.GetFile(), including the Transfer Syntax
    writer.SetImage( reader2.GetImage() );

    gdcm::DataSet &ds = reader1.GetFile().GetDataSet();

    // Make sure that SOPInstanceUID are different
    // Simply removing it is sufficient as gdcm::ImageWriter will generate one by default
    // if not found.
    ds.Remove( gdcm::Tag(0x0008,0x0018) );
    if( !writer.Write() )
    {
        return 1;
    }

    return 0;
}

```

## 29.105 MetaImageMD5Activiz.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;
using gdcm;

/*
 * $ export MONO_PATH=/usr/lib/cli/Activiz.NET:/usr/lib/cli/Kitware.mummy.Runtime-1.0
 * $ mono ./bin/MetaImageMD5Activiz.exe gdcmData/012345.002.050.dcm
 */
public class MetaImageMD5Activiz
{
    public static int ProcessOneMHDMD5(string filename)
    {
        vtkGDCMImageReader reader = vtkGDCMImageReader.
            New();
        reader.FileLowerLeftOn();
        reader.DebugOff();
        int canread = reader.CanReadFile( filename );
        if( canread == 0 )
        {
            string refms = gdcm.Testing.GetMediaStorageFromFile(filename);
            if( gdcm.MediaStorage.IsImage( gdcm.
                MediaStorage.GetMSType(refms) ) )
            {
                System.Console.Write( "Problem with file: " + filename + "\n" );
                return 1;
            }
            // not an image
            return 0;
        }

        reader.SetFileName( filename );
        reader.Update();

        // System.Console.Write(reader.GetOutput());

        vtkMetaImageWriter writer = vtkMetaImageWriter.New();
        writer.SetCompression( false );
        writer.SetInput( reader.GetOutput() );
        string subdir = "MetaImageMD5Activiz";
        string tmpdir = gdcm.Testing.GetTempDirectory( subdir );
        if( !gdcm.PosixEmulation.FileIsDirectory( tmpdir ) )
        {
            gdcm.PosixEmulation.MakeDirectory( tmpdir );
        }
        string mhdfile = gdcm.Testing.GetTempFilename( filename, subdir );

        string rawfile = mhdfile;
        mhdfile += ".mhd";
        rawfile += ".raw";
        writer.SetFileName( mhdfile );
        writer.Write();

        string digestmhd = gdcm.Testing.ComputeFileMD5( mhdfile );
        string digestraw = gdcm.Testing.ComputeFileMD5( rawfile );

        string mhdref = vtkGDCMTesting.GetMHDMD5FromFile(filename);
        string rawref = vtkGDCMTesting.GetRAWMD5FromFile(filename);

        if( mhdref != digestmhd )
        {
            System.Console.Write( "Problem with mhd file: " + filename + "\n" );
            System.Console.Write( digestmhd );
            System.Console.Write( "\n" );
            System.Console.Write( mhdref );
            System.Console.Write( "\n" );
        }
    }
}

```

```

        return 1;
    }
    if( rawref != digestraw )
    {
        System.Console.WriteLine( "Problem with raw file: " + filename + "\n" );
        System.Console.WriteLine( digestraw );
        System.Console.WriteLine( "\n" );
        System.Console.WriteLine( rawref );
        System.Console.WriteLine( "\n" );
        return 1;
    }

    return 0;
}
public static int Main(string[] args)
{
    if ( args.Length == 1 )
    {
        string filename = args[0];
        return ProcessOneMHDMD5( filename );
    }

    // Loop over all gdcmData
    gdcm.Trace.DebugOff();
    gdcm.Trace.WarningOff();
    gdcm.Trace.ErrorOff();

    uint n = gdcm.Testing.GetNumberOfFileNames();
    int ret = 0;
    for( uint i = 0; i < n; ++i )
    {
        string filename = gdcm.Testing.GetFileName( i );
        ret += ProcessOneMHDMD5( filename );
    }
    return ret;
}
}

```

## 29.106 MIPViewer.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
import vtk.*;
import gdcm.*;
import java.io.File;
import java.awt.Canvas;

/*
 * Compilation:
 * CLASSPATH=vtkgdcm.jar:/usr/share/java/vtk.jar javac MIPViewer.java
 *
 * Usage:
 * LD_LIBRARY_PATH=/usr/lib/jvm/java-6-openjdk/jre/lib/amd64/xawt:/usr/lib/jni:. CLASSPATH=/usr/share/java/vtk.jar:vtkgdcm.jar:gdcm.jar:. java MIPViewer BRAINX
 */
public class MIPViewer extends Canvas
{
    static {
        // VTK
        System.loadLibrary("vtkCommonJava");
        System.loadLibrary("vtkFilteringJava");
        System.loadLibrary("vtkIOJava");
        System.loadLibrary("vtkImagingJava");
        System.loadLibrary("vtkGraphicsJava");
        System.loadLibrary("vtkRenderingJava");
        System.loadLibrary("vtkVolumeRenderingJava"); // vtkSmartVolumeMapper
    }
}

```

```

    System.loadLibrary("vtkWidgetsJava"); // vtkBoxWidget
    // VTK-GDCM
    System.loadLibrary("vtkgdcmJava");
}

static FilenamesType fns = new FilenamesType();

protected native int Lock();

protected native int UnLock();

public static void process(String path)
{
    fns.add( path );
}

// Process only files under dir
public static void visitAllFiles(File dir)
{
    if (dir.isDirectory())
    {
        String[] children = dir.list();
        for (int i=0; i<children.length; i++)
        {
            visitAllFiles(new File(dir, children[i]));
        }
    }
    else
    {
        process(dir.getPath());
    }
}

public static void main(String[] args) throws Exception
{
    String dirname = args[0];
    if( !PosixEmulation.FileIsDirectory( dirname ) )
    {
        return;
    }

    File dir = new File(dirname);
    visitAllFiles(dir);

    IPPSorter ipp = new IPPSorter();
    ipp.SetComputeZSpacing( true );
    ipp.SetZSpacingTolerance( 1e-3 );
    boolean b = ipp.Sort( fns );
    if(!b)
    {
        throw new Exception("Could not scan");
    }
    double ippzspacing = ipp.GetZSpacing();

    FilenamesType sorted = ipp.GetFilenames();
    vtkStringArray files = new vtkStringArray();
    long nfiles = sorted.size();
    //for( String f : sorted )
    for (int i = 0; i < nfiles; i++) {
        String f = sorted.get(i);
        files.InsertNextValue( f );
    }
    vtkGDCMImageReader reader = new vtkGDCMImageReader();
    reader.SetFileNames( files );
    reader.Update(); // get spacing value

    double[] spacing = reader.GetOutput().GetSpacing();

    vtkImageChangeInformation change = new vtkImageChangeInformation();
    change.SetInputConnection( reader.GetOutputPort() );
    change.SetOutputSpacing( spacing[0], spacing[1], ippzspacing );

    // Create our volume and mapper
    vtkVolume volume = new vtkVolume();
    vtkSmartVolumeMapper mapper = new vtkSmartVolumeMapper();

    vtkRenderWindowInteractor iren = new vtkRenderWindowInteractor();

    // Add a box widget if the clip option was selected
    vtkBoxWidget box = new vtkBoxWidget();
    box.SetInteractor(iren);

```



```

box.SetPlaceFactor(1.01);
box.SetInput(change.GetOutput());

//box.SetDefaultRenderer(renderer);
box.InsideOutOn();
box.PlaceWidget();
//vtkBoxWidgetCallback callback = vtkBoxWidgetCallback::New();
//callback.SetMapper(mapper);
//box.AddObserver(vtkCommand::InteractionEvent, callback);
//callback.Delete();
// Lock();
// box.EnabledOn();
// Unlock();
box.GetSelectedFaceProperty().SetOpacity(0.0);

mapper.SetInputConnection( change.GetOutputPort() );

// Create our transfer function
vtkColorTransferFunction colorFun = new vtkColorTransferFunction();
vtkPiecewiseFunction opacityFun = new vtkPiecewiseFunction();

// Create the property and attach the transfer functions
vtkVolumeProperty property = new vtkVolumeProperty();
property.IndependentComponentsOn();
property.SetColor( colorFun );
property.SetScalarOpacity( opacityFun );
property.SetInterpolationTypeToLinear();

// connect up the volume to the property and the mapper
volume.SetProperty( property );
volume.SetMapper( mapper );

vtkMedicalImageProperties medprop = reader.GetMedicalImageProperties();
int n = medprop.GetNumberOfWindowLevelPresets();
double opacityWindow = 4096;
double opacityLevel = 2048;

// Override default with value from DICOM files:
for( int i = 0; i < n; ++i )
{
    double wl[] = medprop.GetNthWindowLevelPreset(i);
    //System.out.println( "W/L: " + wl[0] + " " + wl[1] );
    opacityWindow = wl[0];
    opacityLevel = wl[1];
}

colorFun.AddRGBSegment(0.0, 1.0, 1.0, 1.0, 255.0, 1.0, 1.0, 1.0 );
opacityFun.AddSegment( opacityLevel - 0.5*opacityWindow, 0.0,
    opacityLevel + 0.5*opacityWindow, 1.0 );
mapper.SetBlendModeToMaximumIntensity();

// Create the RenderWindow, Renderer
vtkRenderer ren1 = new vtkRenderer();
vtkRenderWindow renWin = new vtkRenderWindow();
renWin.AddRenderer(ren1);

// Set the default window size
renWin.SetSize(600,600);

// Add the volume to the scene
ren1.AddVolume( volume );
ren1.ResetCamera();

iren.SetRenderWindow( renWin );

// interact with data
renWin.Render();

iren.Start();
}
}

```

## 29.107 MpegVideoInfo.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

```

Copyright (c) 2006-2011 Mathieu Malaterre  
 All rights reserved.  
 See Copyright.txt or <http://gdc.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even  
 the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR  
 PURPOSE. See the above copyright notice for more information.

```
=====*/
/*
 * This examples takes in a MPEG2 and write out a Video Endoscopic Image Storage
 * encoded using MPEG2 @ Main Profile
 * ref: http://chrisa.wordpress.com/2007/11/21/decoding-mpeg2-information/
 * See also:
 * http://dvd.sourceforge.net/dvdinfo/mpeghdrs.html#gop
 * http://cvs.linux.hr/cgi-bin/viewcvs.cgi/mpeg\_mod/README.informpeg?view=markup
 * http://www.guru-group.fi/~too/sw/m2vmp2cut/mpeg2info.c
 */

/*
 * Provides information about an MPEG2 file, including the duration, frame rate, aspect
 * ratio, and resolution. Good information about the MPEG2 file structure that helps
 * explain parts of the code can be found here:
 * http://dvd.sourceforge.net/dvdinfo/mpeghdrs.html#gop
 *
 * Copyright (c) 2007 Chris Anderson (chrisa@wordpress.com)
 *
 * This library is free software; you can redistribute it and/or
 * modify it under the terms of the GNU Lesser General Public
 * License as published by the Free Software Foundation; either
 * version 2 of the License, or (at your option) any later version.
 *
 * This library is distributed in the hope that it will be useful,
 * but WITHOUT ANY WARRANTY; without even the implied warranty of
 * MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the GNU
 * Lesser General Public License for more details.
 */
using System;
using System.IO;
using gdc;

public class Mpeg2VideoInfo
{
    #region Member Variables
    private TimeSpan m_startTime = TimeSpan.Zero;
    private TimeSpan m_endTime = TimeSpan.Zero;
    private TimeSpan m_duration = TimeSpan.Zero;
    private eAspectRatios m_aspectRatio = eAspectRatios.Invalid;
    private eFrameRates m_frameRate = 0;
    private int m_pictureWidth = 0;
    private int m_pictureHeight = 0;
    #endregion

    #region Constants
    private const byte PADDING_PACKET = 0xBE;
    private const byte VIDEO_PACKET = 0xE0;
    private const byte AUDIO_PACKET = 0xC0;
    private const byte SYSTEM_PACKET = 0xBB;
    private const byte TIMESTAMP_PACKET = 0xB8;
    private const byte HEADER_PACKET = 0xB3;

    private const int BUFFER_SIZE = 8162; // 8K buffer

    private readonly static TimeSpan EMPTY_TIMESPAN = new TimeSpan(0, 0, -1);
    #endregion

    #region Enumerations
    public enum eFrameRates
    {
        Invalid,
        PulldownNTSC, // 24000d/1001d = 23.976 Hz
        Film, // 24 Hz
        PAL, // 25 Hz
        NTSC, // 30000d/1001d = 29.97 Hz
        DropFrameNTSC, // 30 Hz
        DoubleRatePAL, // 50 Hz
        DoubleRateNTSC, // 59.97 Hz
        DoubleRateDropFrameNTSC // 60 Hz
    }
}
```

```

public enum eAspectRatios
{
    Invalid,
    VGA,           // 1/1
    StandardTV,    // 4/3
    LargeTV,       // 16/9
    Cinema         // 2.21/1
}
#endregion

#region Constructor
public Mpeg2VideoInfo(string file)
{
    ParseMpeg(file);
}
#endregion

#region Public Properties
public TimeSpan StartTime
{
    get { return m_startTime; }
}

public TimeSpan EndTime
{
    get { return m_endTime; }
}

public TimeSpan Duration
{
    get { return m_duration; }
}

public eAspectRatios AspectRatio
{
    get { return m_aspectRatio; }
}

public eFrameRates FrameRate
{
    get { return m_frameRate; }
}

public int PictureWidth
{
    get { return m_pictureWidth; }
}

public int PictureHeight
{
    get { return m_pictureHeight; }
}
#endregion

#region Private Functions
private void ParseMpeg(string file)
{
    FileStream fs = new FileStream(file, FileMode.Open, FileAccess.Read, FileShare.ReadWrite);
    BinaryReader br = new BinaryReader(fs);

    m_startTime = GetStartTimeStampInfo(br);
    m_endTime = GetEndTimeStampInfo(br);

    m_duration = m_endTime.Subtract(m_startTime);

    GetHeaderInfo(br);

    br.Close();
    fs.Close();
}

private TimeSpan GetStartTimeStampInfo(BinaryReader br)
{
    TimeSpan startTime = EMPTY_TIMESPAN;
    byte[] buffer = new byte[BUFFER_SIZE];

    br.BaseStream.Seek(0, SeekOrigin.Begin);

    while (startTime == EMPTY_TIMESPAN && br.BaseStream.Position < br.BaseStream.Length)
    {
        int bytesRead = br.Read(buffer, 0, BUFFER_SIZE);
    }
}

```

```

        for (int offset = 0; offset < readBytes - 8; offset++)
        {
            if (IsStreamMarker(ref buffer, offset, TIMESTAMP_PACKET))
            {
                offset += 4; // Move to the data position which follows the stream header
                uint timeStampEncoded = GetData(ref buffer, offset);
                startTime = DecodeTimeStamp(timeStampEncoded);

                if (startTime != EMPTY_TIMESPAN)
                    break;
            }
        }
    }
    return startTime;
}

private TimeSpan GetEndTimeStampInfo(BinaryReader br)
{
    TimeSpan endTime = EMPTY_TIMESPAN;
    byte[] buffer = new byte[BUFFER_SIZE];

    br.BaseStream.Seek(-BUFFER_SIZE, SeekOrigin.End);

    while (endTime == EMPTY_TIMESPAN && br.BaseStream.Position > BUFFER_SIZE)
    {
        int readBytes = br.Read(buffer, 0, BUFFER_SIZE);

        for (int offset = readBytes - 8; offset >= 0; offset--)
        {
            if (IsStreamMarker(ref buffer, offset, TIMESTAMP_PACKET))
            {
                offset += 4; // Move to the data position which follows the stream header
                uint timeStampEncoded = GetData(ref buffer, offset);
                endTime = DecodeTimeStamp(timeStampEncoded);

                if (endTime != EMPTY_TIMESPAN)
                    break;
            }
        }

        br.BaseStream.Seek(-BUFFER_SIZE * 2, SeekOrigin.Current);
    }

    return endTime;
}

private TimeSpan DecodeTimeStamp(uint timeStampEncoded)
{
    TimeSpan timeStamp = EMPTY_TIMESPAN;

    // Mask out the bits containing the property we are after, then
    // shift the data to the right to get its value
    int hour = (int)(timeStampEncoded & 0x7C000000) >> 26; // Bits 31 -> 27
    int minute = (int)(timeStampEncoded & 0x03F00000) >> 20; // Bits 26 -> 21
    int second = (int)(timeStampEncoded & 0x0007E000) >> 13; // Bits 19 -> 14
    int frame = (int)(timeStampEncoded & 0x00001F80) >> 7; // Bits 13 -> 8 - not used, but included
    for completeness

    timeStamp = new TimeSpan(hour, minute, second);
    return timeStamp;
}

private void GetHeaderInfo(BinaryReader br)
{
    byte[] buffer = new byte[BUFFER_SIZE];

    br.BaseStream.Seek(0, SeekOrigin.Begin);
    br.Read(buffer, 0, BUFFER_SIZE);

    for (int offset = 0; offset < buffer.Length - 4; offset++)
    {
        if (IsStreamMarker(ref buffer, offset, HEADER_PACKET))
        {
            offset += 4; // Move to the data position which follows the stream header
            uint headerData = GetData(ref buffer, offset);

            // Mask out the bits containing the property we are after, then
            // shift the data to the right to get its value
            m_pictureWidth = (int)(headerData & 0xFFFF0000) >> 20;
        }
    }
}

```

```

        m_pictureHeight = (int)(headerData & 0x000FFF00) >> 8;

        uint aspectRatioIndex = (headerData & 0x000000F0) >> 4;
        uint fpsIndex = headerData & 0x0000000F;

        m_aspectRatio = (eAspectRatios)fpsIndex;
        m_frameRate = (eFrameRates)fpsIndex;

        break;
    }
}

private uint GetData(ref byte[] buffer, int offset)
{
    return (uint) ((buffer[offset] << 24) |
        (buffer[offset + 1] << 16) |
        (buffer[offset + 2] << 8) |
        (buffer[offset + 3]));
}

private bool IsStreamMarker(ref byte[] buffer, int offset, byte markerType)
{
    return (buffer[offset] == 0x00 &&
        buffer[offset + 1] == 0x00 &&
        buffer[offset + 2] == 0x01 &&
        buffer[offset + 3] == markerType);
}
#endregion
public static int Main(string[] args)
{
    string file1 = args[0];
    Mpeg2VideoInfo info = new Mpeg2VideoInfo(file1);
    System.Console.WriteLine( info.StartTime );
    System.Console.WriteLine( info.EndTime );
    System.Console.WriteLine( info.Duration );
    System.Console.WriteLine( info.AspectRatio );
    System.Console.WriteLine( info.FrameRate );
    System.Console.WriteLine( info.PictureWidth );
    System.Console.WriteLine( info.PictureHeight );

    ImageReader r = new ImageReader();
    //Image image = new Image();
    Image image = r.GetImage();
    image.SetNumberOfDimensions( 3 );
    DataElement pixeldata = new DataElement( new gdcm.Tag(0x7fe0,0x0010) );

    System.IO.FileStream infile =
        new System.IO.FileStream(file1, System.IO.FileMode.Open, System.IO.FileAccess.Read);
    uint fsize = gdcm.PosixEmulation.FileSize(file1);

    byte[] jstream = new byte[fsize];
    infile.Read(jstream, 0, jstream.Length);

    SmartPtrFrag sq = SequenceOfFragments.New();
    Fragment frag = new Fragment();
    frag.SetByteValue( jstream, new gdcm.VL( (uint)jstream.Length) );
    sq.AddFragment( frag );
    pixeldata.SetValue( sq.__ref__() );

    // insert:
    image.SetDataElement( pixeldata );

    PhotometricInterpretation pi = new PhotometricInterpretation( PhotometricInterpretation.PIType.
        YBR_PARTIAL_420 );
    image.SetPhotometricInterpretation( pi );
    // FIXME hardcoded:
    PixelFormat pixeltype = new PixelFormat(3,8,8,7);
    image.SetPixelFormat( pixeltype );

    // FIXME hardcoded:
    TransferSyntax ts = new TransferSyntax( TransferSyntax.TSType.MPEG2MainProfile);
    image.SetTransferSyntax( ts );

    image.SetDimension(0, (uint)info.PictureWidth);
    image.SetDimension(1, (uint)info.PictureHeight);
    image.SetDimension(2, 721);

    ImageWriter writer = new ImageWriter();
    gdcm.File file = writer.GetFile();
    file.GetHeader().SetDataSetTransferSyntax( ts );

```

```

Anonymizer anon = new Anonymizer();
anon.SetFile( file );

MediaStorage ms = new MediaStorage( MediaStorage.MSType.VideoEndoscopicImageStorage);

UIDGenerator gen = new UIDGenerator();
anon.Replace( new Tag(0x0008,0x16), ms.GetString() );
anon.Replace( new Tag(0x0018,0x40), "25" );
anon.Replace( new Tag(0x0018,0x1063), "40.000000" );
anon.Replace( new Tag(0x0028,0x34), "4\\3" );
anon.Replace( new Tag(0x0028,0x2110), "01" );

writer.SetImage( image );
writer.SetFileName( "dummy.dcm" );
if( !writer.Write() )
{
    System.Console.WriteLine( "Could not write" );
    return 1;
}

return 0;
}
}

```

## 29.108 MPRViewer.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
import vtk.*;
import gdcm.*;
import java.io.File;

/*
 * Compilation:
 * CLASSPATH=vtkgdcm.jar:/usr/share/java/vtk.jar javac MPRViewer.java
 *
 * Usage:
 * LD_LIBRARY_PATH=/usr/lib/jvm/java-6-openjdk/jre/lib/amd64/xawt:/usr/lib/jni:. CLASSPATH=/usr/share/java/
   vtk.jar:vtkgdcm.jar:gdcm.jar:. java MPRViewer BRAINX
 *
 */
public class MPRViewer
{
    static {
        // VTK
        System.loadLibrary("vtkCommonJava");
        System.loadLibrary("vtkFilteringJava");
        System.loadLibrary("vtkIOJava");
        System.loadLibrary("vtkImagingJava");
        System.loadLibrary("vtkGraphicsJava");
        System.loadLibrary("vtkRenderingJava");
        // VTK-GDCM
        System.loadLibrary("vtkgdcmJava");
    }

    static FilenamesType fns = new FilenamesType();

    public static void process(String path)
    {
        fns.add( path );
    }

    // Process only files under dir
    public static void visitAllFiles(File dir)
    {
        if (dir.isDirectory())

```

```

        {
            String[] children = dir.list();
            for (int i=0; i<children.length; i++)
            {
                visitAllFiles(new File(dir, children[i]));
            }
        }
    else
    {
        process(dir.getPath());
    }
}

public static void main(String[] args) throws Exception
{
    String dirname = args[0];
    if( !PosixEmulation.FileIsDirectory( dirname ) )
    {
        return;
    }

    File dir = new File(dirname);
    visitAllFiles(dir);

    IPPSorter ipp = new IPPSorter();
    ipp.SetComputeZSpacing( true );
    ipp.SetZSpacingTolerance( 1e-3 );
    boolean b = ipp.Sort( fns );
    if(!b)
    {
        throw new Exception("Could not scan");
    }
    double ippzspacing = ipp.GetZSpacing();

    FilenamesType sorted = ipp.GetFilenames();
    vtkStringArray files = new vtkStringArray();
    long nfiles = sorted.size();
    //for( String f : sorted )
    for (int i = 0; i < nfiles; i++) {
        String f = sorted.get(i);
        files.InsertNextValue( f );
    }
    vtkGDCMImageReader reader = new vtkGDCMImageReader();
    reader.SetFileNames( files );
    reader.Update(); // get spacing value

    double[] spacing = reader.GetOutput().GetSpacing();

    vtkImageChangeInformation change = new vtkImageChangeInformation();
    change.SetInputConnection( reader.GetOutputPort() );
    change.SetOutputSpacing( spacing[0], spacing[1], ippzspacing );

    // A simple vtkInteractorStyleImage example for
    // 3D image viewing with the vtkImageResliceMapper.
    //
    // Drag Left mouse button to window/level
    // Shift-Left drag to rotate (oblique slice)
    // Shift-Middle drag to slice through image
    // OR Ctrl-Right drag to slice through image

    // Create the RenderWindow, Renderer
    vtkRenderer ren1 = new vtkRenderer();
    vtkRenderWindow renWin = new vtkRenderWindow();
    renWin.AddRenderer(ren1);

    vtkImageResliceMapper im = new vtkImageResliceMapper();
    im.SetInputConnection(change.GetOutputPort());
    im.SliceFacesCameraOn();
    im.SliceAtFocalPointOn();
    im.BorderOff();

    vtkImageProperty ip = new vtkImageProperty();
    ip.SetColorWindow(2000);
    ip.SetColorLevel(1000);
    ip.SetAmbient(0.0);
    ip.SetDiffuse(1.0);
    ip.SetOpacity(1.0);
    ip.SetInterpolationTypeToLinear();

    vtkImageSlice ia = new vtkImageSlice();
    ia.SetMapper(im);

```

```

    ia.SetProperty(ip);

    renl.AddViewProp(ia);
    renl.SetBackground(0.1,0.2,0.4);
    renWin.SetSize(300,300);

    vtkRenderWindowInteractor iren = new vtkRenderWindowInteractor();
    vtkInteractorStyleImage style = new vtkInteractorStyleImage();
    style.SetInteractionModeToImage3D();
    iren.SetInteractorStyle(style);
    renWin.SetInteractor(iren);

    // render the image
    renWin.Render();
    vtkCamera caml = renl.GetActiveCamera();
    caml.ParallelProjectionOn();
    renl.ResetCameraClippingRange();
    renWin.Render();

    iren.Start();
}
}

```

## 29.109 MPRViewer2.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
import vtk.*;
import gdcml.*;
import java.io.File;

/*
 * Compilation:
 * CLASSPATH=vtkgdcml.jar:/usr/share/java/vtk.jar javac MPRViewer2.java
 *
 * Usage:
 * LD_LIBRARY_PATH=/usr/lib/jvm/java-6-openjdk/jre/lib/amd64/xawt:/usr/lib/jni:. CLASSPATH=/usr/share/java/
   vtk.jar:vtkgdcml.jar:gdcml.jar:. java MPRViewer2 BRAINX
 *
 */
public class MPRViewer2
{
    static {
        // VTK
        System.loadLibrary("vtkCommonJava");
        System.loadLibrary("vtkFilteringJava");
        System.loadLibrary("vtkIOJava");
        System.loadLibrary("vtkImagingJava");
        System.loadLibrary("vtkGraphicsJava");
        System.loadLibrary("vtkRenderingJava");
        System.loadLibrary("vtkHybridJava");
        System.loadLibrary("vtkWidgetsJava");
        // VTK-GDCM
        System.loadLibrary("vtkgdcmlJava");
    }

    static FilenamesType fns = new FilenamesType();

    public static void process(String path)
    {
        fns.add( path );
    }

    // Process only files under dir
    public static void visitAllFiles(File dir)
    {

```



```

    if (dir.isDirectory())
    {
        String[] children = dir.list();
        for (int i=0; i<children.length; i++)
        {
            visitAllFiles(new File(dir, children[i]));
        }
    }
    else
    {
        process(dir.getPath());
    }
}

public void dointer(vtkImagePlaneWidget current_widget)
{
    int cstat = current_widget.GetCursorDataStatus();
    double[] v = current_widget.GetCurrentCursorPosition();
    //System.out.println( cstat );
    //System.out.println( v[0] );
    //System.out.println( v[1] );
    //System.out.println( v[2] );
    planeWidgetX.SetSliceIndex( (int)v[0] );
    planeWidgetY.SetSliceIndex( (int)v[1] );
    planeWidgetZ.SetSliceIndex( (int)v[2] );
    planeWidgetX.GetCurrentRenderer().ResetCameraClippingRange();
    planeWidgetY.GetCurrentRenderer().ResetCameraClippingRange();
    planeWidgetZ.GetCurrentRenderer().ResetCameraClippingRange();
}

public void startinterX()
{
    dointer( planeWidgetX );
}

public void interX()
{
    dointer( planeWidgetX );
}

public void endinterX()
{
}

public void startinterY()
{
    dointer( planeWidgetY );
}

public void interY()
{
    dointer( planeWidgetY );
}

public void endinterY()
{
}

public void startinterZ()
{
    dointer( planeWidgetZ );
}

public void interZ()
{
    dointer( planeWidgetZ );
}

public void endinterZ()
{
    //System.out.println( "endinter" );
}

public static void AlignCamera(int slice_number, vtkImagePlaneWidget current_widget)
{
    vtkImageData image = (vtkImageData)current_widget.GetInput();
    vtkRenderer ren = current_widget.GetCurrentRenderer();
    double[] origin = image.GetOrigin();
    double ox = origin[0];
    double oy = origin[1];
    double oz = origin[2];

    int wextent[] = image.GetWholeExtent();
    int xmin = wextent[0];
    int xmax = wextent[1];
    int ymin = wextent[2];
    int ymax = wextent[3];
    int zmin = wextent[4];
    int zmax = wextent[5];
}

```

```

double[] spacing = image.GetSpacing();
double sx = spacing[0];
double sy = spacing[1];
double sz = spacing[2];

double cx = ox+(0.5*(xMax-xMin))*sx;
double cy = oy+(0.5*(yMax-yMin))*sy;
double cz = oz+(0.5*(zMax-zMin))*sz;
double vx = 0, vy = 0, vz = 0;
double nx = 0, ny = 0, nz = 0;
int iaxis = current_widget.GetPlaneOrientation();
if ( iaxis == 0 ) {
    vz = -1;
    nx = ox + xMax*sx;
    cx = ox + slice_number*sx;
}
else if ( iaxis == 1 ) {
    vz = -1;
    ny = oy+yMax*sy;
    cy = oy+slice_number*sy;
}
else {
    vy = 1;
    nz = oz+zMax*sz;
    cz = oz+slice_number*sz;
}
double px = cx+nx*2;
double py = cy+ny*2;
double pz = cz+nz*3;

vtkCamera camera = ren.GetActiveCamera();
camera.SetViewUp(vx, vy, vz);
camera.SetFocalPoint(cx, cy, cz);
camera.SetPosition(px, py, pz);
camera.OrthogonalizeViewUp();
ren.ResetCameraClippingRange();
}

private vtkImagePlaneWidget planeWidgetX = new vtkImagePlaneWidget();
private vtkImagePlaneWidget planeWidgetY = new vtkImagePlaneWidget();
private vtkImagePlaneWidget planeWidgetZ = new vtkImagePlaneWidget();

public void config()
{
    //System.out.println( "config" );
    planeWidgetX.GetCurrentRenderer().ResetCamera();
    planeWidgetY.GetCurrentRenderer().ResetCamera();
    planeWidgetZ.GetCurrentRenderer().ResetCamera();
}

public void Run(String dirname)
{
    File dir = new File(dirname);
    visitAllFiles(dir);

    IPPSorter ipp = new IPPSorter();
    ipp.SetComputeZSpacing( true );
    ipp.SetZSpacingTolerance( 1e-3 );
    boolean b = ipp.Sort( fns );
    if(!b)
    {
        //throw new Exception("Could not scan");
    }
    double ippzspacing = ipp.GetZSpacing();

    FilenamesType sorted = ipp.GetFilenames();
    vtkStringArray files = new vtkStringArray();
    long nfiles = sorted.size();
    //for( String f : sorted )
    for (int i = 0; i < nfiles; i++) {
        String f = sorted.get(i);
        files.InsertNextValue( f );
    }
    vtkGDCMImageReader reader = new vtkGDCMImageReader();
    reader.SetFileNames( files );
    reader.Update(); // get spacing value

    double[] spacing = reader.GetOutput().GetSpacing();

    vtkImageChangeInformation change = new vtkImageChangeInformation();
    change.SetInputConnection( reader.GetOutputPort() );

```

```
change.SetOutputSpacing( spacing[0], spacing[1], ippszspacing );
change.Update();

System.out.println( change.GetOutput().toString() );

vtkRenderer ren1 = new vtkRenderer();
ren1.SetViewport(0., 0., 0.333, 1);
ren1.SetBackground(0.1,0.2,0.4);
vtkRenderer ren2 = new vtkRenderer();
ren2.SetViewport(0.333, 0., 0.667, 1);
ren2.SetBackground(0.1,0.2,0.4);
vtkRenderer ren3 = new vtkRenderer();
ren3.SetViewport(0.667, 0., 1., 1.);
ren3.SetBackground(0.1,0.2,0.4);

vtkRenderWindow renWin = new vtkRenderWindow();
renWin.AddRenderer(ren1);
renWin.AddRenderer(ren2);
renWin.AddRenderer(ren3);

vtkRenderWindowInteractor iren = new vtkRenderWindowInteractor();
iren.SetRenderWindow(renWin);

vtkInteractorStyleImage style = new vtkInteractorStyleImage();
iren.SetInteractorStyle( style );

vtkCellPicker picker = new vtkCellPicker();
picker.SetTolerance(0.005);

vtkProperty ipwProp = new vtkProperty();

//vtkImagePlaneWidget planeWidgetX = new vtkImagePlaneWidget();
planeWidgetX.SetInteractor(iren);
planeWidgetX.SetCurrentRenderer(ren1);
planeWidgetX.SetDefaultRenderer(ren1);
planeWidgetX.RestrictPlaneToVolumeOn();
planeWidgetX.SetTexturePlaneProperty(ipwProp);
//planeWidgetX.GetPlaneProperty().SetColor(1,0,0);
//planeWidgetX.TextureInterpolateOff();
//planeWidgetX.SetResliceInterpolateToNearestNeighbour();
planeWidgetX.SetInput( change.GetOutput() );
planeWidgetX.SetPlaneOrientationToXAxes();
planeWidgetX.SetSliceIndex(62);
planeWidgetX.SetPicker(picker);
planeWidgetX.SetKeyPressActivationValue('x');
planeWidgetX.On();
planeWidgetX.InteractionOn();

//vtkImagePlaneWidget planeWidgetY = new vtkImagePlaneWidget();
planeWidgetY.SetInteractor(iren);
planeWidgetY.SetCurrentRenderer(ren2);
planeWidgetY.SetDefaultRenderer(ren2);
planeWidgetY.RestrictPlaneToVolumeOn();
planeWidgetY.SetTexturePlaneProperty(ipwProp);
//planeWidgetY.GetPlaneProperty().SetColor(1,0,0);
//planeWidgetY.TextureInterpolateOff();
//planeWidgetY.SetResliceInterpolateToNearestNeighbour();
planeWidgetY.SetInput( change.GetOutput() );
planeWidgetY.SetLookupTable( planeWidgetX.GetLookupTable() );
planeWidgetY.SetPlaneOrientationToYAxes();
planeWidgetY.SetSliceIndex(32);
planeWidgetY.SetPicker(picker);
planeWidgetY.SetKeyPressActivationValue('y');
planeWidgetY.On();

//vtkImagePlaneWidget planeWidgetZ = new vtkImagePlaneWidget();
planeWidgetZ.SetInteractor(iren);
planeWidgetZ.SetCurrentRenderer(ren3);
planeWidgetZ.SetDefaultRenderer(ren3);
planeWidgetZ.RestrictPlaneToVolumeOn();
planeWidgetZ.SetTexturePlaneProperty(ipwProp);
//planeWidgetZ.GetPlaneProperty().SetColor(1,0,0);
//planeWidgetZ.TextureInterpolateOff();
//planeWidgetZ.SetResliceInterpolateToNearestNeighbour();
planeWidgetZ.SetInput( change.GetOutput() );
planeWidgetZ.SetLookupTable( planeWidgetX.GetLookupTable() );
planeWidgetZ.SetPlaneOrientationToZAxes();
planeWidgetZ.SetSliceIndex(32);
planeWidgetZ.SetPicker(picker);
planeWidgetZ.SetKeyPressActivationValue('z');
```

```

planeWidgetZ.On();

iren.Initialize();

renWin.Render();
AlignCamera(52, planeWidgetX);
AlignCamera(32, planeWidgetY);
AlignCamera(32, planeWidgetZ);

planeWidgetX.GetCurrentRenderer().ResetCamera();
planeWidgetY.GetCurrentRenderer().ResetCamera();
planeWidgetZ.GetCurrentRenderer().ResetCamera();

renWin.Render();

planeWidgetX.AddObserver("StartInteractionEvent", this,"startinterX");
planeWidgetX.AddObserver("InteractionEvent", this,"interX");
planeWidgetX.AddObserver("EndInteractionEvent", this,"endinterX");
planeWidgetY.AddObserver("StartInteractionEvent", this,"startinterY");
planeWidgetY.AddObserver("InteractionEvent", this,"interY");
planeWidgetY.AddObserver("EndInteractionEvent", this,"endinterY");
planeWidgetZ.AddObserver("StartInteractionEvent", this,"startinterZ");
planeWidgetZ.AddObserver("InteractionEvent", this,"interZ");
planeWidgetZ.AddObserver("EndInteractionEvent", this,"endinterZ");

iren.AddObserver("ConfigureEvent", this,"config");

iren.Start();
}

public static void main(String[] args) throws Exception
{
    String dirname = args[0];
    if( !PosixEmulation.FileIsDirectory( dirname ) )
    {
        return;
    }

    MPRViewer2 me = new MPRViewer2();
    me.Run( dirname );
}
}

```

## 29.110 MrProtocol.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
*
*/

/*
28 - 'MrProtocol' VM 1, VR UN, SyngoDT 0, NoOfItems 6, Data '### ASCCONV BEGIN ###
ulVersion                               = 0xbee332
tSequenceFileName                       = "%SiemensSeq%\fl_fq_shphs"
tProtocolName                           = "flash+AF8-100+AF8-through-plane+AF8-V"
tReferenceImage0                         = "1.3.12.2.1107.5.2.9.16041.30000007062106100181200004658"
tReferenceImage1                         = "1.3.12.2.1107.5.2.9.16041.30000007062106100181200004635"
tReferenceImage2                         = "1.3.12.2.1107.5.2.9.16041.30000007062106100181200004683"
ucScanRegionPosValid                    = 0x1
sProtConsistencyInfo.tBaselineString    = "N4_VB11A_LATEST_20031004"
sProtConsistencyInfo.flNominalB0        = 1.494
sProtConsistencyInfo.flGMax              = 22
sProtConsistencyInfo.flRiseTime          = 10
sGRADSPEC.sEddyCompensationX.aflAmplitude[0] = 0.0141111
sGRADSPEC.sEddyCompensationX.aflAmplitude[1] = 0.057038

```

```
sGRADSPEC.sEddyCompensationX.aflAmplitude[2] = -0.00986504
sGRADSPEC.sEddyCompensationX.aflAmplitude[3] = 0.00247627
sGRADSPEC.sEddyCompensationX.aflAmplitude[4] = 0.0026377
sGRADSPEC.sEddyCompensationX.aflTimeConstant[0] = 1.53826
sGRADSPEC.sEddyCompensationX.aflTimeConstant[1] = 0.746617
sGRADSPEC.sEddyCompensationX.aflTimeConstant[2] = 0.339236
sGRADSPEC.sEddyCompensationX.aflTimeConstant[3] = 0.0309809
sGRADSPEC.sEddyCompensationX.aflTimeConstant[4] = 0.00067694
sGRADSPEC.sEddyCompensationY.aflAmplitude[0] = 0.0156411
sGRADSPEC.sEddyCompensationY.aflAmplitude[1] = 0.0440623
sGRADSPEC.sEddyCompensationY.aflAmplitude[2] = -0.00782663
sGRADSPEC.sEddyCompensationY.aflAmplitude[3] = 0.00186828
sGRADSPEC.sEddyCompensationY.aflAmplitude[4] = 0.00154504
sGRADSPEC.sEddyCompensationY.aflTimeConstant[0] = 1.47145
sGRADSPEC.sEddyCompensationY.aflTimeConstant[1] = 0.750538
sGRADSPEC.sEddyCompensationY.aflTimeConstant[2] = 0.339397
sGRADSPEC.sEddyCompensationY.aflTimeConstant[3] = 0.0312962
sGRADSPEC.sEddyCompensationY.aflTimeConstant[4] = 0.000895133
sGRADSPEC.sEddyCompensationZ.aflAmplitude[0] = 0.00618504
sGRADSPEC.sEddyCompensationZ.aflAmplitude[1] = 0.00313121
sGRADSPEC.sEddyCompensationZ.aflAmplitude[2] = 0.000289346
sGRADSPEC.sEddyCompensationZ.aflAmplitude[3] = -0.00019677
sGRADSPEC.sEddyCompensationZ.aflAmplitude[4] = 7.66445e-005
sGRADSPEC.sEddyCompensationZ.aflTimeConstant[0] = 3.37462
sGRADSPEC.sEddyCompensationZ.aflTimeConstant[1] = 0.999351
sGRADSPEC.sEddyCompensationZ.aflTimeConstant[2] = 0.0174646
sGRADSPEC.sEddyCompensationZ.aflTimeConstant[3] = 0.0110094
sGRADSPEC.sEddyCompensationZ.aflTimeConstant[4] = 0.00199922
sGRADSPEC.bEddyCompensationValid = 1
sGRADSPEC.sB0CompensationX.aflAmplitude[0] = 0.307474
sGRADSPEC.sB0CompensationX.aflAmplitude[1] = 0.029337
sGRADSPEC.sB0CompensationX.aflAmplitude[2] = -0.187118
sGRADSPEC.sB0CompensationX.aflTimeConstant[0] = 0.98583
sGRADSPEC.sB0CompensationX.aflTimeConstant[1] = 0.0308443
sGRADSPEC.sB0CompensationX.aflTimeConstant[2] = 0.000466792
sGRADSPEC.sB0CompensationY.aflAmplitude[0] = 0.365257
sGRADSPEC.sB0CompensationY.aflAmplitude[1] = -0.318647
sGRADSPEC.sB0CompensationY.aflAmplitude[2] = -0.0118978
sGRADSPEC.sB0CompensationY.aflTimeConstant[0] = 0.61535
sGRADSPEC.sB0CompensationY.aflTimeConstant[1] = 0.488831
sGRADSPEC.sB0CompensationY.aflTimeConstant[2] = 0.00199991
sGRADSPEC.sB0CompensationZ.aflAmplitude[0] = -0.44647
sGRADSPEC.sB0CompensationZ.aflAmplitude[1] = -0.0455154
sGRADSPEC.sB0CompensationZ.aflAmplitude[2] = -0.0304901
sGRADSPEC.sB0CompensationZ.aflTimeConstant[0] = 0.959231
sGRADSPEC.sB0CompensationZ.aflTimeConstant[1] = 0.0720189
sGRADSPEC.sB0CompensationZ.aflTimeConstant[2] = 0.00190141
sGRADSPEC.bB0CompensationValid = 1
sGRADSPEC.sCrossTermCompensationXY.aflAmplitude[0] = 0.00105046
sGRADSPEC.sCrossTermCompensationXY.aflTimeConstant[0] = 0.842014
sGRADSPEC.sCrossTermCompensationXZ.aflAmplitude[0] = -0.00150189
sGRADSPEC.sCrossTermCompensationXZ.aflTimeConstant[0] = 0.736169
sGRADSPEC.sCrossTermCompensationYX.aflAmplitude[0] = -5.5278e-005
sGRADSPEC.sCrossTermCompensationYX.aflTimeConstant[0] = 0.228697
sGRADSPEC.sCrossTermCompensationYZ.aflAmplitude[0] = 0.000307999
sGRADSPEC.sCrossTermCompensationYZ.aflTimeConstant[0] = 1.19431
sGRADSPEC.sCrossTermCompensationZX.aflAmplitude[0] = -0.000286868
sGRADSPEC.sCrossTermCompensationZX.aflTimeConstant[0] = 0.665979
sGRADSPEC.sCrossTermCompensationZY.aflAmplitude[0] = 0.000355175
sGRADSPEC.sCrossTermCompensationZY.aflTimeConstant[0] = 0.844189
sGRADSPEC.bCrossTermCompensationValid = 1
sGRADSPEC.lOffsetX = 25
sGRADSPEC.lOffsetY = 84
sGRADSPEC.lOffsetZ = 47
sGRADSPEC.bOffsetValid = 1
sGRADSPEC.lDelayX = 12
sGRADSPEC.lDelayY = 11
sGRADSPEC.lDelayZ = 9
sGRADSPEC.bDelayValid = 1
sGRADSPEC.flSensitivityX = 0.000264087
sGRADSPEC.flSensitivityY = 0.000272009
sGRADSPEC.flSensitivityZ = 0.000272677
sGRADSPEC.bSensitivityValid = 1
sGRADSPEC.alShimCurrent[0] = 183
sGRADSPEC.alShimCurrent[1] = -25
sGRADSPEC.alShimCurrent[2] = -85
sGRADSPEC.alShimCurrent[3] = 378
sGRADSPEC.alShimCurrent[4] = 82
sGRADSPEC.bShimCurrentValid = 1
sGRADSPEC.ucMode = 0x2
sTXSPEC.asNucleusInfo[0].tNucleus = "1H"
```

```

sTXSPEC.asNucleusInfo[0].lFrequency      = 63684693
sTXSPEC.asNucleusInfo[0].bFrequencyValid = 1
sTXSPEC.asNucleusInfo[0].flReferenceAmplitude = 359.734
sTXSPEC.asNucleusInfo[0].bReferenceAmplitudeValid = 1
sTXSPEC.asNucleusInfo[0].flAmplitudeCorrection = 1
sTXSPEC.asNucleusInfo[0].bAmplitudeCorrectionValid = 1
sTXSPEC.asNucleusInfo[1].bFrequencyValid = 1
sTXSPEC.asNucleusInfo[1].bReferenceAmplitudeValid = 1
sTXSPEC.asNucleusInfo[1].bAmplitudeCorrectionValid = 1
sTXSPEC.arFPULSE[0].tName                 = "03GreFCE"
sTXSPEC.arFPULSE[0].bAmplitudeValid       = 0x1
sTXSPEC.arFPULSE[0].flAmplitude           = 147.095
sTXSPEC.arFPULSE[1].tName                 = "02GreFCE"
sTXSPEC.arFPULSE[1].bAmplitudeValid       = 0x1
sTXSPEC.arFPULSE[1].flAmplitude           = 147.095
sTXSPEC.arFPULSE[2].tName                 = "01GreFCE"
sTXSPEC.arFPULSE[2].bAmplitudeValid       = 0x1
sTXSPEC.arFPULSE[2].flAmplitude           = 147.095
sTXSPEC.lNoOfTraPulses                    = 3
sTXSPEC.lBTB1ParallelCapacity              = 2
sTXSPEC.lBTB1SerialCapacity                = 24
sTXSPEC.lBTB2ParallelCapacity              = 2
sTXSPEC.lBTB2SerialCapacity                = 26
sTXSPEC.bBTBValid                          = 1
sTXSPEC.flKDynMagnitudeMin                 = 0.5
sTXSPEC.flKDynMagnitudeMax                 = 1.5
sTXSPEC.flKDynMagnitudeClipLow             = 0.96
sTXSPEC.flKDynMagnitudeClipHigh            = 1.04
sTXSPEC.flKDynPhaseMax                     = 0.698132
sTXSPEC.flKDynPhaseClip                    = 0.174533
sTXSPEC.bKDynValid                         = 1
sTXSPEC.ucRFPulseType                      = 0x1
sTXSPEC.ucExcitMode                        = 0x1
sTXSPEC.ucSimultaneousExcitation            = 0x1
sRXSPEC.lGain                              = 1
sRXSPEC.bGainValid                         = 1
sRXSPEC.aFFT_SCALE[0].lRxChannel            = 1
sRXSPEC.aFFT_SCALE[0].flFactor              = 1.06857
sRXSPEC.aFFT_SCALE[0].bValid                = 1
sRXSPEC.aFFT_SCALE[1].lRxChannel            = 2
sRXSPEC.aFFT_SCALE[1].flFactor              = 1.07454
sRXSPEC.aFFT_SCALE[1].bValid                = 1
sRXSPEC.aFFT_SCALE[2].lRxChannel            = 3
sRXSPEC.aFFT_SCALE[2].flFactor              = 1.06622
sRXSPEC.aFFT_SCALE[2].bValid                = 1
sRXSPEC.aFFT_SCALE[3].lRxChannel            = 4
sRXSPEC.aFFT_SCALE[3].flFactor              = 1.06524
sRXSPEC.aFFT_SCALE[3].bValid                = 1
sRXSPEC.aFFT_SCALE[4].lRxChannel            = 5
sRXSPEC.aFFT_SCALE[4].flFactor              = 0.982692
sRXSPEC.aFFT_SCALE[4].bValid                = 1
sRXSPEC.aFFT_SCALE[5].lRxChannel            = 6
sRXSPEC.aFFT_SCALE[5].flFactor              = 0.988603
sRXSPEC.aFFT_SCALE[5].bValid                = 1
sRXSPEC.aFFT_SCALE[6].lRxChannel            = 7
sRXSPEC.aFFT_SCALE[6].flFactor              = 0.981538
sRXSPEC.aFFT_SCALE[6].bValid                = 1
sRXSPEC.aFFT_SCALE[7].lRxChannel            = 8
sRXSPEC.aFFT_SCALE[7].flFactor              = 1.00856
sRXSPEC.aFFT_SCALE[7].bValid                = 1
sRXSPEC.bVariCapVoltagesValid              = 1
sRXSPEC.alDwellTime[0]                     = 8500
sAdjFreSpec.ulMode                          = 0x1
sAdjFreSpec.ucAdjWithBC                     = 0x1
sAdjTraSpec.ucAdjWithBC                     = 0x1
sAdjShimSpec.ulMode                         = 0x1
sAdjShimSpec.ucAdjWithBC                     = 0x1
sAdjWatSupSpec.ulMode                       = 0x1
sAdjWatSupSpec.ucAdjWithBC                  = 0x1
alTR[0]                                     = 37000
lContrasts                                 = 1
alTE[0]                                    = 4000
acFlowComp[0]                              = 1
lCombinedEchoes                             = 1
sSliceArray.asSlice[0].sPosition.dSag       = 35.31199581
sSliceArray.asSlice[0].sPosition.dCor       = -8.387765754
sSliceArray.asSlice[0].sPosition.dTra       = -23.13178296
sSliceArray.asSlice[0].sNormal.dSag         = 0.771051253
sSliceArray.asSlice[0].sNormal.dCor         = 0.5863890019
sSliceArray.asSlice[0].sNormal.dTra         = -0.2482496801
sSliceArray.asSlice[0].dThickness           = 6

```

```

sSliceArray.asSlice[0].dPhaseFOV      = 187.5
sSliceArray.asSlice[0].dReadoutFOV    = 250
sSliceArray.lSize                      = 1
sSliceArray.lSag                      = 1
sSliceArray.lConc                     = 1
sSliceArray.ucMode                    = 0x1
sSliceArray.sTSat.dThickness          = 40
sSliceArray.sTSat.dGap                = 10
sGroupArray.asGroup[0].nSize          = 1
sGroupArray.asGroup[0].dDistFact      = 0.2
sGroupArray.anMember[1]               = -1
sGroupArray.lSize                     = 1
sGroupArray.sPSat.dThickness          = 50
sGroupArray.sPSat.dGap                = 10
sAutoAlign.dAAMatrix[0]               = 1
sAutoAlign.dAAMatrix[5]               = 1
sAutoAlign.dAAMatrix[10]              = 1
sAutoAlign.dAAMatrix[15]              = 1
sNavigatorPara.ucRespComp              = 0x4
sPrepPulses.ucFatSat                  = 0x4
sPrepPulses.ucWaterSat                = 0x4
sPrepPulses.ucInversion                = 0x4
sPrepPulses.ucSatRecovery              = 0x1
sPrepPulses.ucFatSatMode               = 0x2
sKSpace.lBaseResolution                = 256
sKSpace.lPhaseEncodingLines            = 192
sKSpace.dPhaseResolution                = 1
sKSpace.lPartitions                    = 32
sKSpace.lImagesPerSlab                 = 32
sKSpace.dSliceResolution                = 1
sKSpace.ucPhasePartialFourier           = 0x10
sKSpace.ucSlicePartialFourier           = 0x10
sKSpace.ucAveragingMode                 = 0x2
sKSpace.ucMultiSliceMode                = 0x1
sKSpace.ucDimension                     = 0x2
sKSpace.ucAsymmetricEchoAllowed         = 0x1
sKSpace.unReordering                   = 0x1
sFastImaging.lEPIFactor                 = 1
sFastImaging.lTurboFactor               = 1
sFastImaging.lSegments                  = 3
sFastImaging.ulEnableRFSpoiling         = 0x1
sPhysioImaging.lSignal1                 = 2
sPhysioImaging.lMethod1                 = 2
sPhysioImaging.lSignal2                 = 1
sPhysioImaging.lMethod2                 = 1
sPhysioImaging.lPhases                  = 21
sPhysioImaging.lRetroGatedImages        = 16
sPhysioImaging.sPhysioECG.lScanWindow   = 805
sPhysioImaging.sPhysioECG.lTriggerPulses = 1
sPhysioImaging.sPhysioECG.lTriggerWindow = 5
sPhysioImaging.sPhysioECG.lArrhythmiaDetection = 1
sPhysioImaging.sPhysioECG.lCardiacGateOnThreshold = 100000
sPhysioImaging.sPhysioECG.lCardiacGateOffThreshold = 700000
sPhysioImaging.sPhysioPulse.lTriggerPulses = 1
sPhysioImaging.sPhysioPulse.lTriggerWindow = 5
sPhysioImaging.sPhysioPulse.lCardiacGateOnThreshold = 100000
sPhysioImaging.sPhysioPulse.lCardiacGateOffThreshold = 700000
sPhysioImaging.sPhysioExt.lTriggerPulses = 1
sPhysioImaging.sPhysioExt.lTriggerWindow = 5
sPhysioImaging.sPhysioExt.lCardiacGateOnThreshold = 100000
sPhysioImaging.sPhysioExt.lCardiacGateOffThreshold = 700000
sPhysioImaging.sPhysioResp.lRespGateThreshold = 20
sPhysioImaging.sPhysioResp.lRespGatePhase = 2
sPhysioImaging.sPhysioResp.dGatingRatio = 0.3
sSpecPara.lPhaseCyclingType             = 1
sSpecPara.lPhaseEncodingType            = 1
sSpecPara.lRFExcitationBandwidth         = 1
sSpecPara.ucRemoveOversampling           = 0x1
sSpecPara.lDecouplingType                = 1
sSpecPara.lNOEType                       = 1
sSpecPara.lExcitationType                 = 1
sSpecPara.lSpectralSuppression            = 1
sDiffusion.ulMode                        = 0x1
sAngio.sFlowArray.asElm[0].nVelocity     = 100
sAngio.sFlowArray.asElm[0].nDir          = 0x4
sAngio.sFlowArray.lSize                  = 1
sAngio.ucPCFlowMode                      = 0x2
sAngio.ucTOFInflow                       = 0x4
sAngio.ucRephasedImage                   = 0x1
sAngio.ucPhaseImage                      = 0x1
sEllipticalFilter.ucMode                  = 0x1

```

```

sPat.lAccelFactPE = 1
sPat.lAccelFact3D = 1
sPat.ucPATMode = 0x1
sPat.ucRefScanMode = 0x1
ucAutoMovie = 0x1
ucDisableChangeStoreImages = 0x1
ucReconstructionMode = 0x1
ucPHAPSMode = 0x1
ucDixon = 0x1
lAverages = 2
adFlipAngleDegree[0] = 30
lScanTimeSec = 103
lTotalScanTimeSec = 112
dRefSNR = 165404.1473
dRefSNR_VOI = 165404.1473
tdefaultEVAProt = "%SiemensEvaDefProt%\Inline\Inline.evp"
tcurrentEVAProt = "%CURRENTEVAPROT%\EVA2A5.tmp"
sCOIL_SELECT_MEAS.asList[0].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[0].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[0].sCoilElementID.tElement = "PP6"
sCOIL_SELECT_MEAS.asList[0].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[0].lRxChannelConnected = 1
sCOIL_SELECT_MEAS.asList[1].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[1].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[1].sCoilElementID.tElement = "PP5"
sCOIL_SELECT_MEAS.asList[1].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[1].lRxChannelConnected = 1
sCOIL_SELECT_MEAS.asList[2].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[2].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[2].sCoilElementID.tElement = "PP3"
sCOIL_SELECT_MEAS.asList[2].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[2].lRxChannelConnected = 2
sCOIL_SELECT_MEAS.asList[3].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[3].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[3].sCoilElementID.tElement = "PP4"
sCOIL_SELECT_MEAS.asList[3].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[3].lRxChannelConnected = 3
sCOIL_SELECT_MEAS.asList[4].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[4].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[4].sCoilElementID.tElement = "PP2"
sCOIL_SELECT_MEAS.asList[4].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[4].lRxChannelConnected = 4
sCOIL_SELECT_MEAS.asList[5].sCoilElementID.tCoilID = "6_Ch_Body_P"
sCOIL_SELECT_MEAS.asList[5].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[5].sCoilElementID.tElement = "PP1"
sCOIL_SELECT_MEAS.asList[5].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[5].lRxChannelConnected = 4
sCOIL_SELECT_MEAS.asList[6].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[6].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[6].sCoilElementID.tElement = "PA6"
sCOIL_SELECT_MEAS.asList[6].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[6].lRxChannelConnected = 5
sCOIL_SELECT_MEAS.asList[7].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[7].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[7].sCoilElementID.tElement = "PA5"
sCOIL_SELECT_MEAS.asList[7].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[7].lRxChannelConnected = 5
sCOIL_SELECT_MEAS.asList[8].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[8].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[8].sCoilElementID.tElement = "PA3"
sCOIL_SELECT_MEAS.asList[8].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[8].lRxChannelConnected = 6
sCOIL_SELECT_MEAS.asList[9].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[9].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[9].sCoilElementID.tElement = "PA4"
sCOIL_SELECT_MEAS.asList[9].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[9].lRxChannelConnected = 7
sCOIL_SELECT_MEAS.asList[10].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[10].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[10].sCoilElementID.tElement = "PA2"
sCOIL_SELECT_MEAS.asList[10].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[10].lRxChannelConnected = 8
sCOIL_SELECT_MEAS.asList[11].sCoilElementID.tCoilID = "6_Ch_Body_A"
sCOIL_SELECT_MEAS.asList[11].sCoilElementID.lCoilCopy = 1
sCOIL_SELECT_MEAS.asList[11].sCoilElementID.tElement = "PA1"
sCOIL_SELECT_MEAS.asList[11].lElementSelected = 1
sCOIL_SELECT_MEAS.asList[11].lRxChannelConnected = 8
sCOIL_SELECT_MEAS.sCOILPLUGS.aulPlugId[0] = 0xff
sCOIL_SELECT_MEAS.sCOILPLUGS.aulPlugId[1] = 0x76
sCOIL_SELECT_MEAS.sCOILPLUGS.aulPlugId[2] = 0x78
sCOIL_SELECT_MEAS.sCOILPLUGS.aulPlugId[3] = 0x87

```



```

sCOIL_SELECT_MEAS.sCOILPLUGS.auiPlugId[4] = 0x67
sCOIL_SELECT_MEAS.sCOILPLUGS.auiNmbrOfNibbles[0] = 0x2
sCOIL_SELECT_MEAS.sCOILPLUGS.auiNmbrOfNibbles[1] = 0x2
sCOIL_SELECT_MEAS.sCOILPLUGS.auiNmbrOfNibbles[2] = 0x2
sCOIL_SELECT_MEAS.sCOILPLUGS.auiNmbrOfNibbles[3] = 0x2
sCOIL_SELECT_MEAS.sCOILPLUGS.auiNmbrOfNibbles[4] = 0x2
sEFISPEC.bEFIDataValid = 1
### ASCCONV END ###
'
*/

/*
 * Table of equivalence:
 *
ulVersion = 0xb332
<=>
27 - 'MrProtocolVersion' VM 1, VR IS, SyngoDT 6, NoOfItems 6, Data '12510002'
*/

#include "gdcmReader.h"
#include "gdcmImageReader.h"
#include "gdcmImageWriter.h"
#include "gdcmCSAHeader.h"
#include "gdcmAttribute.h"
#include "gdcmGlobal.h"
#include "gdcmDicts.h"

#include <map>

#include <math.h>

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcm::ImageReader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }

    gdcm::CSAHeader csa;
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    //const gdcm::PrivateTag &t1 = csa.GetCSAImageHeaderInfoTag();
    const gdcm::PrivateTag &t2 = csa.GetCSASeriesHeaderInfoTag();

    if( ds.FindDataElement( t2 ) )
    {
        csa.LoadFromDataElement( ds.GetDataElement( t2 ) );
        //csa.Print( std::cout );
    }

    if( !csa.FindCSAElementByName( "MrProtocol" ) )
    {
        return 1;
    }
    const gdcm::CSAElement &csael = csa.GetCSAElementByName( "MrProtocol" );
    //std::cout << csael << std::endl;

    const gdcm::ByteValue *bv = csael.GetByteValue();
    if( !bv )
    {
        return 1;
    }
    std::string str(bv->GetPointer(), bv->GetLength());
    std::istringstream is(str);
    std::string s;
    typedef std::map< std::string, std::string > MyMapType;
    MyMapType mymap;
    while( std::getline(is, s) )
    {
        std::string::size_type pos = s.find( '=' );
        if( pos != std::string::npos )
        {
            std::string sub1 = s.substr(0, pos);
            sub1.erase( sub1.find_last_not_of( ' ' ) + 1 );
            std::string sub2 = s.substr(pos+1); // skip the '=' char

```

```

        sub2.erase( 0, sub2.find_first_not_of(' '));
        //std::cout << sub1 << std::endl;
        mymap.insert( MyMapType::value_type(sub1, sub2) );
    }
    else
    {
        // ### ASCCONV BEGIN ###
        // ### ASCCONV END ###
    }
}

const char fourierstr[] = "sKSpace.ucSlicePartialFourier";
const gdcmm::CSAHeaderDict &csadict =
    gdcmm::Global::GetInstance().GetDicts().
    GetCSAHeaderDict();
const gdcmm::CSAHeaderDictEntry &fourier = csadict.
    GetCSAHeaderDictEntry( fourierstr );
std::cout << fourier << std::endl;
MyMapType::const_iterator it = mymap.find ( fourierstr );
if( it == mymap.end() ) return 1;
//std::cout << it->second << std::endl;
const std::string &partial_fourier = it->second;
if( partial_fourier == "0x1" )
{
    std::cout << "partial fourier is 4/8" << std::endl;
}
else if( partial_fourier == "0x2" )
{
    std::cout << "partial fourier is 5/8" << std::endl;
}
else if( partial_fourier == "0x4" )
{
    std::cout << "partial fourier is 6/8" << std::endl;
}
else if( partial_fourier == "0x8" )
{
    std::cout << "partial fourier is 7/8" << std::endl;
}
else if( partial_fourier == "0x10" )
{
    std::cout << "partial fourier is 8/8" << std::endl;
}
else
{
    std::cerr << "Impossible: " << partial_fourier << std::endl;
    return 1;
}
}

/*
This is the Flip Angle:
adFlipAngleDegree[0]                = 30

One can find it also in the protocol:

...
    <ParamFunctor."<TlmapFunctor">">
    {
        <Class> "<TlmapFunctor@IceImagePostProcFunctors">

        <ParamBool."<EXECUTE">"> { }
        <ParamDouble."<Flip_deg">"> { <Precision> 16 14.7378520000000000 }
    }
...

*/
// Below is an attempt to play with the CSAHeader dict:
#if 0
const char gspec[] = "sGRADSPEC.flSensitivityX";
it = mymap.find( gspec );
if( it == mymap.end() ) return 1;
const std::string &dummy = it->second;
std::cout << dummy << std::endl;

const gdcmm::CSAHeaderDictEntry &csaentry = csadict.
    GetCSAHeaderDictEntry( gspec );
std::cout << csaentry << std::endl;
#endif

/*
sSliceArray.ucMode -- should be in (1, 2, 4)
enum SeriesMode
{
    ASCENDING    = 0x01,

```

```

        DESCENDING = 0x02,
        INTERLEAVED = 0x04
    };
*/
const char sliceorderstr[] = "sSliceArray.ucMode";
const gdcm::CSAHeaderDictEntry &sliceorder = csadict.
    GetCSAHeaderDictEntry( sliceorderstr );
std::cout << sliceorder << std::endl;

it = mymap.find ( sliceorderstr );
if( it == mymap.end() ) return 1;
const std::string &slice_order = it->second;
if( slice_order == "0x1" )
{
    std::cout << "slice_order: ASCENDING" << std::endl;
}
else if( slice_order == "0x2" )
{
    std::cout << "slice_order: DESCENDING" << std::endl;
}
else if( slice_order == "0x4" )
{
    std::cout << "slice_order: INTERLEAVED" << std::endl;
}
else
{
    std::cerr << "Impossible: " << slice_order << std::endl;
    return 1;
}

return 0;
}

```

## 29.111 NewSequence.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/NewSequence.exe gdcmData/012345.002.050.dcm out.dcm
 */
using System;
//using gdcm;

public class NewSequence
{
    public static byte[] StrToByteArray(string str)
    {
        System.Text.ASCIIEncoding encoding=new System.Text.ASCIIEncoding();
        return encoding.GetBytes(str);
    }

    public static int Main(string[] argv)
    {
        string file1 = argv[0];
        string file2 = argv[1];

        gdcm.Reader r = new gdcm.Reader();
        r.SetFileName( file1 );
        if ( ! r.Read() )
        {
            return 1;
        }
    }
}

```

```

gdcM.File f = r.GetFile();
gdcM.DataSet ds = f.GetDataSet();
// tsis = gdcM.Tag(0x0008,0x2112) # SourceImageSequence

// Create a dataelement
gdcM.DataElement de = new gdcM.DataElement(new
    gdcM.Tag(0x0010, 0x2180));
string occ = "Occupation";
de.SetByteValue( StrToByteArray(occ), new gdcM.VL((uint)occ.Length));
de.SetVR(new gdcM.VR(gdcM.VR.VRType.SH));

// Create an item
gdcM.Item it = new gdcM.Item();
it.SetVLToUndefined(); // Needed to not popup error message
//it.InsertDataElement(de)
gdcM.DataSet nds = it.GetNestedDataSet();
nds.Insert(de);

// Create a Sequence
gdcM.SmartPtrSQ sq = gdcM.SequenceOfItems.New();
sq.SetLengthToUndefined();
sq.AddItem(it);

// Insert sequence into data set
gdcM.DataElement des = new gdcM.DataElement(new
    gdcM.Tag(0x0400,0x0550));
des.SetVR(new gdcM.VR(gdcM.VR.VRType.SQ));
des.SetValue(sq.__ref__());
des.SetVLToUndefined();

ds.Insert(des);

gdcM.Writer w = new gdcM.Writer();
w.SetFile( f );
w.SetFileName( file2 );
if ( !w.Write() )
    return 1;

return 0;
}

```

## 29.112 NewSequence.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python NewSequence.py input.dcm output.dcm
19
20
21 Thanks to Robert Irie for code
22 """
23
24 import sys
25 import gdcM
26
27 if __name__ == "__main__":
28
29     file1 = sys.argv[1]
30     file2 = sys.argv[2]
31
32     r = gdcM.Reader()
33     r.SetFileName( file1 )

```

```

34  if not r.Read():
35      sys.exit(1)
36
37  f = r.GetFile()
38  ds = f.GetDataSet()
39  #tsis = gdcM.Tag(0x0008,0x2112) # SourceImageSequence
40
41  # Create a dataelement
42  de = gdcM.DataElement(gdcM.Tag(0x0010, 0x2180))
43  de.SetByteValue("Occupation", gdcM.VL(len("Occupation")))
44  de.SetVR(gdcM.VR(gdcM.VR.SH))
45
46  # Create an item
47  it=gdcM.Item()
48  it.SetVLToUndefined() # Needed to not popup error message
49  #it.InsertDataElement(de)
50  nds=it.GetNestedDataSet()
51  nds.Insert(de)
52
53  # Create a Sequence
54  sq=gdcM.SequenceOfItems().New()
55  sq.SetLengthToUndefined()
56  sq.AddItem(it)
57
58  # Insert sequence into data set
59  des=gdcM.DataElement(gdcM.Tag(0x0400,0x0550))
60  des.SetVR(gdcM.VR(gdcM.VR.SQ))
61  des.SetValue(sq.__ref__())
62  des.SetVLToUndefined()
63
64  ds.Insert(des)
65
66  w = gdcM.Writer()
67  w.SetFile( f )
68  w.SetFileName( file2 )
69  if not w.Write():
70      sys.exit(1)

```

## 29.113 offscreenimage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMImageReader.h"
#include "vtkRenderWindow.h"
#include "vtkRenderer.h"
#include "vtkImageMapToWindowLevelColors.h"
#include "vtkImageActor.h"
#include "vtkPNGWriter.h"
#include "vtkWindowToImageFilter.h"
#include "vtkMedicalImageProperties.h"

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        return 1;
    }
    const char *filename = argv[1];

    vtkGDCMImageReader *reader = vtkGDCMImageReader::New();
    reader->SetFileName( filename );
    reader->Update(); // important to read the window/level info

    vtkMedicalImageProperties *prop = reader->GetMedicalImageProperties();

    vtkRenderWindow *renWin = vtkRenderWindow::New();

```

```

renWin->OffScreenRenderingOn();

vtkRenderer *renderer = vtkRenderer::New();
renWin->AddRenderer(renderer);

vtkImageMapToWindowLevelColors *windowlevel = vtkImageMapToWindowLevelColors::New();
#if (VTK_MAJOR_VERSION >= 6)
windowlevel->SetInputConnection( reader->GetOutputPort() );
#else
windowlevel->SetInput( reader->GetOutput() );
#endif
unsigned int n = prop->GetNumberOfWindowLevelPresets();
if( n )
{
    // Take the first one by default:
    const double *wl = prop->GetNthWindowLevelPreset(0);
    windowlevel->SetWindow( wl[0] );
    windowlevel->SetLevel( wl[1] );
}

vtkImageActor *actor = vtkImageActor::New();
#if (VTK_MAJOR_VERSION >= 6)
actor->SetInputData( windowlevel->GetOutput() );
#else
actor->SetInput( windowlevel->GetOutput() );
#endif

renderer->AddActor( actor );

renWin->Render();

vtkWindowToImageFilter *w2if = vtkWindowToImageFilter::New();
w2if->SetInput( renWin );

vtkPNGWriter *wr = vtkPNGWriter::New();
#if (VTK_MAJOR_VERSION >= 6)
wr->SetInputConnection( w2if->GetOutputPort() );
#else
wr->SetInput( w2if->GetOutput() );
#endif
wr->SetFileName( "offscreenimage.png" );
wr->Write();

reader->Delete();
renWin->Delete();
renderer->Delete();
windowlevel->Delete();
actor->Delete();
w2if->Delete();
wr->Delete();

return 0;
}

```

## 29.114 PatchFile.cxx

This is a C++ example on how to use [gdcm::Attribute](#)

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * The image was a broken file where the Pixel Data element was 8 times too big
 * Apparently multiplying the BitsAllocated to 4 and multiplying the number of
 * frames by 2 would solve the problem
 */

```

```

* This C++ code can be used to patch the header.
*/

#include "gdcmReader.h"
#include "gdcmImageReader.h"
#include "gdcmWriter.h"
#include "gdcmDataSet.h"
#include "gdcmAttribute.h"

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        return 1;
    }
    const char *f = argv[1];
    const char *out = argv[2];
    gdcm::Reader r;
    r.SetFileName( f );
    if( !r.Read() )
    {
        return 1;
    }

    gdcm::File &file = r.GetFile();
    gdcm::DataSet& ds = file.GetDataSet();
    // (0028,0100) US 16 # 2, 1 BitsAllocated
    // (0028,0101) US 16 # 2, 1 BitsStored
    // (0028,0102) US 15 # 2, 1 HighBit
    //
    {
        gdcm::Attribute<0x28,0x100> at;
        at.SetFromDataElement( ds.GetDataElement( at.
            GetTag() ) );
        if( at.GetValue() != 8 )
        {
            return 1;
        }
        at.SetValue( 32 );
        ds.Replace( at.GetAsDataElement() );
    }
    {
        gdcm::Attribute<0x28,0x101> at;
        at.SetFromDataElement( ds.GetDataElement( at.
            GetTag() ) );
        if( at.GetValue() != 8 )
        {
            return 1;
        }
        at.SetValue( 32 );
        ds.Replace( at.GetAsDataElement() );
    }
    {
        gdcm::Attribute<0x28,0x102> at;
        at.SetFromDataElement( ds.GetDataElement( at.
            GetTag() ) );
        if( at.GetValue() != 7 )
        {
            return 1;
        }
        at.SetValue( 31 );
        ds.Replace( at.GetAsDataElement() );
    }
    // (0028,0008) IS [56] # 2, 1 NumberOfFrames

    {
        gdcm::Attribute<0x28,0x8> at;
        at.SetFromDataElement( ds.GetDataElement( at.
            GetTag() ) );
        at.SetValue( at.GetValue() * 2 );
        ds.Replace( at.GetAsDataElement() );
    }

    gdcm::Writer w;
    w.SetFile( file );
    w.SetCheckFileMetaInformation( false );
    w.SetFileName( out );
    if( !w.Write() )
    {
        return 1;
    }
}

```

```
// Now let's see if we can read it as an image:
gdcmm::ImageReader ir;
ir.SetFileName( out );
if(!ir.Read())
{
    return 1;
}
gdcmm::Image &image = ir.GetImage();
unsigned long len = image.GetBufferLength();
const gdcmm::ByteValue *bv = ir.GetFile().GetDataSet().
    GetDataElement( gdcmm::Tag(0x7fe0,0x0010) ).GetByteValue();
if( !bv || len != bv->GetLength() )
{
    return 1;
}
std::cout << bv->GetLength() << " " << len << std::endl;

std::cout << "Success to rewrite image !" << std::endl;
image.Print( std::cout );
return 0;
}
```

## 29.115 PhilipsPrivateRescaleInterceptSlope.py

```
1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python
19 """
20
21 import gdcmm
22 import sys
23
24 filename = sys.argv[1]
25 tmpfile = "/tmp/philips_rescaled.dcm"
26
27
28 # Need to access some private tags, read the file :
29 reader = gdcmm.Reader()
30 reader.SetFileName( filename )
31 if not reader.Read():
32     sys.exit(1)
33
34 ds = reader.GetFile().GetDataSet()
35
36 #print ds
37 # (2005,1409) DS 4 0.0
38 # (2005,140a) DS 16 1.52283272283272
39
40 # (2005,0014) LO 26 Philips MR Imaging DD 005
41 tag1 = gdcmm.PrivateTag(0x2005,0x09,"Philips MR Imaging DD 005")
42 tag2 = gdcmm.PrivateTag(0x2005,0x0a,"Philips MR Imaging DD 005")
43 print tag1
44 print tag2
45
46 # make sure to do a copy, we want the private tag to remain
47 # otherwise gdcmm gives us a reference
48 e11 = gdcmm.DataElement( ds.GetDataElement( tag1 ) )
49 print e11
50 e12 = gdcmm.DataElement( ds.GetDataElement( tag2 ) )
51 print e12
52
```



```

53 # (0028,1052) DS [-1000]                # 6, 1 RescaleIntercept
54 # (0028,1053) DS [1]                    # 2, 1 RescaleSlope
55
56 el1.SetTag( gdcm.Tag(0x0028,0x1052) )
57 el2.SetTag( gdcm.Tag(0x0028,0x1053) )
58
59 ds.Insert( el1 )
60 ds.Insert( el2 )
61
62 w = gdcm.Writer()
63 w.SetCheckFileMetaInformation( False )
64 w.SetFileName( tmpfile )
65 w.SetFile( reader.GetFile() )
66 if not w.Write():
67     sys.exit(1)
68
69 print "success"

```

## 29.116 PlaySound.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python PlaySound.py input.dcm
19 """
20
21 import gdcm
22 import sys
23
24 #filename = "/home/mmalaterre/Creatis/gdcmDataExtra/gdcmNonImageData/audio_from_rafael_sanguinetti.dcm"
25 filename = sys.argv[1]
26 print filename
27
28 r = gdcm.Reader()
29 r.SetFileName( filename )
30 if not r.Read():
31     sys.exit(1)
32
33 ds = r.GetFile().GetDataSet()
34
35 waveformtag = gdcm.Tag(0x5400,0x0100)
36 waveformsq = ds.GetDataElement( waveformtag )
37 #print waveformsq
38
39 #print dir(waveformsq)
40
41 items = waveformsq.GetSequenceOfItems()
42
43 if not items.GetNumberOfItems():
44     sys.exit(1)
45
46 item = items.GetItem(1)
47 #print item
48
49 waveformds = item.GetNestedDataSet()
50 #print waveformds
51
52 waveformdatatag = gdcm.Tag(0x5400,0x01010)
53 waveformdata = waveformds.GetDataElement( waveformdatatag )
54
55 #print waveformdata.GetPointer()
56 bv = waveformdata.GetByteValue()
57 print dir(bv)

```

```

58
59 #print bv.GetPointer()
60 print bv.GetLength()
61 l = 116838
62
63 file='test.wav'
64 myfile = open(file, "wb")
65 s = bv.GetPointer()
66 for i in range(0, l):
67     myfile.write(s[i])
68 myfile.close()
69
70 # http://mail.python.org/pipermail/python-list/2004-October/288905.html
71 if sys.platform.startswith('win'):
72     from winsound import PlaySound, SND_FILENAME, SND_ASYNC
73     PlaySound(file, SND_FILENAME|SND_ASYNC)
74 elif sys.platform.find('linux')>-1:
75     from wave import open as waveOpen
76     from ossaudiodev import open as ossOpen
77     s = waveOpen(file,'rb')
78     (nc,sw,fr,nf,comptype, compname) = s.getparams( )
79     dsp = ossOpen('/dev/dsp','w')
80     try:
81         from ossaudiodev import AFMT_S16_NE
82     except ImportError:
83         if byteorder == "little":
84             AFMT_S16_NE = ossaudiodev.AFMT_S16_LE
85         else:
86             AFMT_S16_NE = ossaudiodev.AFMT_S16_BE
87     dsp.setparameters(AFMT_S16_NE, nc, fr)
88     data = s.readframes(nf)
89     s.close()
90     dsp.write(data)
91     dsp.close()

```

## 29.117 pmsct\_rgb1.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example shows how to rewrite a ELSCINT1/PMSCT_RGB1 compressed
 * image so that it is readable by most 3rd party software (DICOM does
 * not specify this particular encoding).
 * This is required for the sake of interoperability with any standard
 * conforming DICOM system.
 *
 * Everything done in this code is for the sole purpose of writing interoperable
 * software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
 * If you believe anything in this code violates any law or any of your rights,
 * please contact us (gdcm-developers@lists.sourceforge.net) so that we can
 * find a solution.
 *
 * Everything you do with this code is at your own risk, since decompression
 * algorithm was not written from specification documents.
 *
 * Special thanks to:
 * Jean-Pierre Roux for providing the sample datasets
 */
#include "gdcmReader.h"
#include "gdcmPrivateTag.h"
#include "gdcmAttribute.h"
#include "gdcmImageWriter.h"

void delta_decode(const unsigned char *data_in, size_t data_size,
                  std::vector<unsigned char> &new_stream, unsigned short pc, size_t w, size_t h)
{

```

```

const size_t plane_size = h * w;
const size_t outputlen = 3 * plane_size;
new_stream.resize( outputlen );

assert( data_size != outputlen );
if( data_size == outputlen )
{
    return;
}
typedef unsigned char byte;
enum {
    COLORMODE   = 0x81,
    ESCMODE     = 0x82,
    REPEATMODE  = 0x83
};

byte* src = (byte*)data_in;
byte* dest = (byte*)&new_stream[0];
union { byte gray; byte rgb[3]; } pixel;
pixel.rgb[0] = pixel.rgb[1] = pixel.rgb[2] = 0;
// always start in grayscale mode
bool graymode = true;
size_t dx = 1;
size_t dy = 3;
// algorithm works with both planar configuration
// It does produce surprising greenish background color for planar
// configuration is 0, while the nested Icon SQ display a nice black
// background
if (pc)
{
    dx = plane_size;
    dy = 1;
}
size_t ps = plane_size;

// The following is highly unoptimized as we have nested if statement in a while loop
// we need to switch from one algorithm to ther other (RGB <-> GRAY)
while (ps)
{
    // next byte:
    byte b = *src++;
    assert( src < data_in + data_size );
    // mode selection:
    switch ( b )
    {
        case ESCMODE:
            // Used to treat a byte 81/82/83 as a normal byte
            if (graymode)
            {
                pixel.gray += *src++;
                dest[0*dx] = pixel.gray;
                dest[1*dx] = pixel.gray;
                dest[2*dx] = pixel.gray;
            }
            else
            {
                pixel.rgb[0] += *src++;
                pixel.rgb[1] += *src++;
                pixel.rgb[2] += *src++;
                dest[0*dx] = pixel.rgb[0];
                dest[1*dx] = pixel.rgb[1];
                dest[2*dx] = pixel.rgb[2];
            }
            dest += dy;
            ps--;
            break;
        case REPEATMODE:
            // repeat mode (RLE)
            b = *src++;
            ps -= b;
            if (graymode)
            {
                while (b-- > 0)
                {
                    dest[0*dx] = pixel.gray;
                    dest[1*dx] = pixel.gray;
                    dest[2*dx] = pixel.gray;
                    dest += dy;
                }
            }
            else

```

```

        {
            while (b-- > 0)
            {
                dest[0*dx] = pixel.rgb[0];
                dest[1*dx] = pixel.rgb[1];
                dest[2*dx] = pixel.rgb[2];
                dest += dy;
            }
        }
        break;
case COLORMODE:
    // We are swithing from one mode to the other. The stream contains an intermixed
    // compression of RGB codec and GRAY codec. Each one not knowing of the other
    // reset old value to 0.
    if (graymode)
    {
        graymode = false;
        pixel.rgb[0] = pixel.rgb[1] = pixel.rgb[2] = 0;
    }
    else
    {
        graymode = true;
        pixel.gray = 0;
    }
    break;
default:
    // This is identical to ESCMODE, it would be nicer to use fall-through
    if (graymode)
    {
        pixel.gray += b;
        dest[0*dx] = pixel.gray;
        dest[1*dx] = pixel.gray;
        dest[2*dx] = pixel.gray;
    }
    else
    {
        pixel.rgb[0] += b;
        pixel.rgb[1] += *src++;
        pixel.rgb[2] += *src++;
        dest[0*dx] = pixel.rgb[0];
        dest[1*dx] = pixel.rgb[1];
        dest[2*dx] = pixel.rgb[2];
    }
    dest += dy;
    ps--;
    break;
    } // end switch
} // end while
}

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
        return 1;
    }
    const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

    // (07a1,1011) CS [PMSCT_RGB1] # 10,1 Tamar Compression Type
    const gdcm::PrivateTag tcompressiontype(0x07a1,0x0011,"ELSCINT1");
    if( !ds.FindDataElement( tcompressiontype ) ) return 1;
    const gdcm::DataElement& compressiontype = ds.GetDataElement(
        tcompressiontype );
    if ( compressiontype.IsEmpty() ) return 1;
    const gdcm::ByteValue *bv = compressiontype.GetByteValue();
    std::string comprle = "PMSCT_RLE1";
    std::string comprgb = "PMSCT_RGB1";
    bool isrle = false;
    bool isrgb = false;
    if( strncmp( bv->GetPointer(), comprle.c_str(), comprle.size() ) == 0 )
    {
        isrle = true;
        return 1;
    }
    if( strncmp( bv->GetPointer(), comprgb.c_str(), comprgb.size() ) == 0 )
    {

```

```

    isrgb = true;
}
if( !isrgb && !isrle ) return 1;

const gdcm::PrivateTag tcompressedpixeldata(0x07a1,0x000a,"ELSCINT1");
if( !ds.FindDataElement( tcompressedpixeldata ) ) return 1;
const gdcm::DataElement& compressionpixeldata = ds.
    GetDataElement( tcompressedpixeldata);
if ( compressionpixeldata.IsEmpty() ) return 1;
const gdcm::ByteValue * bv2 = compressionpixeldata.GetByteValue();

gdcm::Attribute<0x0028,0x0006> at0;
at0.SetFromDataSet( ds );
gdcm::Attribute<0x0028,0x0010> at1;
at1.SetFromDataSet( ds );
gdcm::Attribute<0x0028,0x0011> at2;
at2.SetFromDataSet( ds );

std::vector<unsigned char> buffer;
delta_decode((const unsigned char*)bv2->GetPointer(), bv2->GetLength(), buffer,
    at0.GetValue(), at1.GetValue(), at2.GetValue() );

gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
pixeldata.SetVR( gdcm::VR::OW );
pixeldata.SetByteValue( (char*)&buffer[0], (uint32_t)buffer.size() );
// TODO we should check that decompress byte buffer match the expected size (row*col*...)

// Add the pixel data element
reader.GetFile().GetDataSet().Replace( pixeldata );

reader.GetFile().GetHeader().SetDataSetTransferSyntax(
    gdcm::TransferSyntax::ExplicitVRLittleEndian);
gdcm::Writer writer;
writer.SetFile( reader.GetFile() );

// Cleanup stuff:
// remove the compressed pixel data:
// FIXME: should I remove more private tags ? all of them ?
// oh well this is just an example
// use gdcm::Anonymizer::RemovePrivateTags if needed...
writer.GetFile().GetDataSet().Remove( compressionpixeldata.
    GetTag() );
std::string outfilename;
if (argc > 2)
    outfilename = argv[2];
else
    outfilename = "outrgb.dcm";
writer.SetFileName( outfilename.c_str() );
if( !writer.Write() )
{
    std::cerr << "Failed to write" << std::endl;
    return 1;
}

std::cout << "success !" << std::endl;

return 0;
}

```

## 29.118 PrivateDict.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 """

```

```

17
18 import gdcM
19 import sys,os
20
21 if __name__ == "__main__":
22     #gdcM.Trace.DebugOn()
23     globInst = gdcM.Global.GetInstance()
24     # Try to load Part3.xml file
25     # This file is too big for being accessible directly at runtime.
26     globInst.LoadResourcesFiles()
27
28
29     # Get a private tag from the runtime dicts. LoadResourcesFiles could
30     # have failed but this has no impact on the private dict
31
32     d = globInst.GetDicts()
33     print d.GetDictEntry( gdcM.Tag(0x0029,0x0010) ,"SIEMENS CSA HEADER" )
34     pd = d.GetPrivateDict()
35     print pd.GetDictEntry( gdcM.PrivateTag(0x0029,0x0010,"SIEMENS CSA HEADER") )

```

## 29.119 PublicDict.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcM.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
/*
 * Dummy example to show GDCM Dict(s) API (Part 6) + Collected Private Attributes:
 */

#include "gdcMGlobal.h"
#include "gdcMDicts.h"
#include "gdcMDict.h"
#include "gdcMCSAHeader.h"
#include "gdcMPrivateTag.h"

int main(int , char *[])
{
    const gdcM::Global& g = gdcM::Global::GetInstance(); // sum of all
        knowledge !
    const gdcM::Dicts &dicts = g.GetDicts();
    const gdcM::Dict &pub = dicts.GetPublicDict(); // Part 6

    //std::cout << pub << std::endl;

    // 3 different ways to access the same information

    // 1. From the public dict only:
    gdcM::Tag patient_name(0x10,0x10);
    const gdcM::DictEntry &entry1 = pub.GetDictEntry(patient_name);
    std::cout << entry1 << std::endl;

    // 2. From all dicts:
    const gdcM::DictEntry &entry2 = dicts.GetDictEntry(patient_name);
    std::cout << entry2 << std::endl;

    // 3. This solution is the most flexible solution as you can request using the same
    // API either a public tag or a private tag
    const char *strowner = 0;
    const gdcM::DictEntry &entry3 = dicts.GetDictEntry(patient_name,strowner);
    std::cout << entry3 << std::endl;

    // Private attributes:

    // try with a private tag now:
    const gdcM::PrivateTag &private_tag =
        gdcM::CSAHeader::GetCSAImageHeaderInfoTag();
    //std::cout << private_tag << std::endl;

```

```

const gdcm::DictEntry &entry4 = dicts.GetDictEntry(private_tag,private_tag.
    GetOwner());
std::cout << entry4 << std::endl;

// Let's pretend that private lookup is on 0x10xx elements:
gdcm::PrivateTag dummy = private_tag;
dummy.SetElement( (uint16_t)(0x1000 + dummy.GetElement()) );
const gdcm::DictEntry &entry5 = dicts.GetDictEntry(dummy,dummy.
    GetOwner());
std::cout << entry5 << std::endl;

return 0;
}

```

## 29.120 QIDO-RS.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
#include "gdcmReader.h"
#include "gdcmWriter.h"
#include "gdcmJSON.h"

/*
 * Simple QIDO-RS round-trip to test implementation of gdcm::JSON
 * See Supl66 for details
 */
int main(int argc, char *argv[])
{
    if( argc < 2 ) return 1;
    using namespace gdcm;
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() ) return 1;

    gdcm::JSON json;
    json.PrettyPrintOn();
    std::stringstream ss;
    const gdcm::File & f = reader.GetFile();
    json.Code( f.GetDataSet(), ss);

    std::cout << ss.str() << std::endl;

    gdcm::Writer w;
    gdcm::File & ff = w.GetFile();
    ff.GetHeader().SetDataSetTransferSyntax(
        gdcm::TransferSyntax::ExplicitVRLittleEndian );
    if( !json.Decode(ss, ff.GetDataSet() ) )
    {
        std::cerr << "Could not decode" << std::endl;
        return 1;
    }
    w.SetFileName( "/tmp/debug.dcm" );
    if( !w.Write() ) return 1;

    return 0;
}

```

## 29.121 ReadAndDumpDICOMDIR.cxx

```

/*=====

```

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre

All rights reserved.

See Copyright.txt or <http://gdcm.sourceforge.net/Copyright.html> for details.

This software is distributed WITHOUT ANY WARRANTY; without even  
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR  
PURPOSE. See the above copyright notice for more information.

```

=====*/
/*
 * This example shows how to read and dump a DICOMDIR File
 *
 * Thanks:
 *   Tom Marynowski (lordglub gmail) for contributing this example
 */
#include "gdcmReader.h"
#include "gdcmMediaStorage.h"

typedef std::set<gdcm::DataElement> DataElementSet;
typedef DataElementSet::const_iterator ConstIterator;

int main(int argc, char *argv [])
{
    if( argc < 2 ) return 1;
    const char *filename = argv[1];

    gdcm::Reader reader;
    reader.SetFileName( filename);
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }
    std::stringstream strm;

    gdcm::File &file = reader.GetFile();
    gdcm::DataSet &ds = file.GetDataSet();
    gdcm::FileMetaInformation &fmi = file.GetHeader();

    gdcm::MediaStorage ms;
    ms.SetFromFile(file);
    if( ms != gdcm::MediaStorage::MediaStorageDirectoryStorage
        )
    {
        std::cout << "This file is not a DICOMDIR" << std::endl;
        return 1;
    }

    if (fmi.FindDataElement( gdcm::Tag (0x0002, 0x0002)))
    {
        strm.str("");
        fmi.GetDataElement( gdcm::Tag (0x0002, 0x0002) ).
            GetValue().Print(strm);
    }
    else
    {
        std::cerr << " Media Storage Sop Class UID not present" << std::endl;
    }

    //TODO il faut trimer strm.str() avant la comparaison au cas ou...
    if ("1.2.840.10008.1.3.10"!=strm.str())
    {
        std::cout << "This file is not a DICOMDIR" << std::endl;
        return 1;
    }

    ConstIterator it = ds.GetDES().begin();

    for( ; it != ds.GetDES().end(); ++it)
    {
        if (it->GetTag()==gdcm::Tag (0x0004, 0x1220))
        {
            const gdcm::DataElement &de = (*it);
            // ne pas utiliser GetSequenceOfItems pour extraire les items
            gdcm::SmartPointer<gdcm::SequenceOfItems> sqi =de.
                GetValueAsSQ();
            unsigned int itemused = 1;

```



```

while (itemused<=sqi->GetNumberOfItems())

{
    strm.str("");

    if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0004, 0x1430)))
        sqi->GetItem(itemused).GetDataElement(gdcm::Tag (0x0004, 0x1430)).
GetValue().Print(strm);

    //TODO il faut trimer strm.str() avant la comparaison
    while((strm.str()=="PATIENT")||((strm.str()=="PATIENT ")))
    {
        std::cout << strm.str() << std::endl;
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0010, 0x0010)))
            sqi->GetItem(itemused).GetDataElement(gdcm::Tag (0x0010, 0x0010))
.GetValue().Print(strm);
        std::cout << "PATIENT NAME : " << strm.str() << std::endl;

        //PATIENT ID
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0010, 0x0020)))
            sqi->GetItem(itemused).GetDataElement(gdcm::Tag (0x0010, 0x0020))
.GetValue().Print(strm);
        std::cout << "PATIENT ID : " << strm.str() << std::endl;

        /*ADD TAG TO READ HERE*/
        std::cout << "===== " << std::endl;
        itemused++;
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0004, 0x1430)))
            sqi->GetItem(itemused).GetDataElement(gdcm::Tag (0x0004, 0x1430))
.GetValue().Print(strm);

        //TODO il faut trimer strm.str() avant la comparaison
        while((strm.str()=="STUDY")||((strm.str()=="STUDY ")))
        {
            std::cout << " " << strm.str() << std::endl;
            //UID
            strm.str("");
            if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0020, 0x000d)))
                sqi->GetItem(itemused).GetDataElement(
gdcm::Tag (0x0020, 0x000d)).GetValue().Print(strm);
            std::cout << " " << strm.str() << std::endl;

            //STUDY DATE
            strm.str("");
            if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0008, 0x0020)))
                sqi->GetItem(itemused).GetDataElement(
gdcm::Tag (0x0008, 0x0020)).GetValue().Print(strm);
            std::cout << " " << strm.str() << std::endl;

            //STUDY DESCRIPTION
            strm.str("");
            if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0008, 0x1030)))
                sqi->GetItem(itemused).GetDataElement(
gdcm::Tag (0x0008, 0x1030)).GetValue().Print(strm);
            std::cout << " " << strm.str() << std::endl;

            /*ADD TAG TO READ HERE*/
            std::cout << " " << "===== " << std::endl;

            itemused++;
            strm.str("");
            if (sqi->GetItem(itemused).FindDataElement(
gdcm::Tag (0x0004, 0x1430)))
                sqi->GetItem(itemused).GetDataElement(
gdcm::Tag (0x0004, 0x1430)).GetValue().Print(strm);

            //TODO il faut trimer strm.str() avant la comparaison
            while((strm.str()=="SERIES")||((strm.str()=="SERIES ")))
            {
                std::cout << " " << strm.str() << std::endl;

```

```

        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0020, 0x000e)))
            sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0020, 0x000e)).GetValue().Print(strm);
        std::cout << "                SERIE UID" << strm.str() << std::endl;

        //SERIE MODALITY
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0008, 0x0060)))
            sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0008, 0x0060)).GetValue().Print(strm);
        std::cout << "                SERIE MODALITY" << strm.str() << std::endl;

        //SERIE DESCRIPTION
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0008, 0x103e)))
            sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0008, 0x103e)).GetValue().Print(strm);
        std::cout << "                SERIE DESCRIPTION" << strm.str() << std::endl;

        /*ADD TAG TO READ HERE*/

        std::cout << "                " << "===== " << std::endl;
        itemused++;
        strm.str("");
        if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0004, 0x1430)))
            sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0004, 0x1430)).GetValue().Print(strm);

        //TODO il faut trimer strm.str() avant la comparaison
        while ((strm.str()=="IMAGE")||((strm.str()=="IMAGE ")))
            // if(tmp=="IMAGE")
            {
                std::cout << "                " << strm.str() << std::endl;

                //UID
                strm.str("");
                if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0004, 0x1511)))
                    sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0004, 0x1511)).GetValue().Print(strm);
                std::cout << "                IMAGE UID : " << strm.str() << std::endl;

                //PATH de l'image
                strm.str("");
                if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0004, 0x1500)))
                    sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0004, 0x1500)).GetValue().Print(strm);
                std::cout << "                IMAGE PATH : " << strm.str() << std::endl;
                /*ADD TAG TO READ HERE*/

                if(itemused < sqi->GetNumberOfItems())
                {itemused++;
                }else{break;}

                strm.str("");

                if (sqi->GetItem(itemused).FindDataElement(
gdcmm::Tag (0x0004, 0x1430)))
                    sqi->GetItem(itemused).GetDataElement (
gdcmm::Tag (0x0004, 0x1430)).GetValue().Print(strm);

            }
        }
    }
    itemused++;
}
}
return 0;

```

```
}

```

## 29.122 ReadAndDumpDICOMDIR.py

```
1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 # File: ReadAndDumpDICOMDIR.py
14 #
15 # Author: Lukas Batteau (lbatteau gmail)
16 #
17 # This example shows how to read and dump a DICOMDIR File.
18 # Based on Tom Marynowski's (lordglub gmail) example.
19 #
20 # Usage:
21 # python ReadAndDumpDICOMDIR.py [DICOMDIR file]
22 #####
23
24
25
26 import sys
27 import gdcm
28
29 if __name__ == "__main__":
30     # Check arguments
31     if (len(sys.argv) < 2):
32         # No filename passed
33         print "No input filename found"
34         quit()
35
36     filename = sys.argv[1]
37
38
39     # Read file
40     reader = gdcm.Reader()
41     reader.SetFileName(filename)
42     if (not reader.Read()):
43         print "Unable to read %s" % (filename)
44         quit()
45
46     file = reader.GetFile()
47
48     # Retrieve header information
49     fileMetaInformation = file.GetHeader()
50     print fileMetaInformation
51
52     # Retrieve data set
53     dataSet = file.GetDataSet()
54     #print dataSet
55
56     # Check media storage
57     mediaStorage = gdcm.MediaStorage()
58     mediaStorage.SetFromFile(file)
59     if (gdcm.MediaStorage.GetMSType(str(mediaStorage)) !=
60 gdcm.MediaStorage.MediaStorageDirectoryStorage):
61         # File is not a DICOMDIR
62         print "This file is not a DICOMDIR (Media storage type: %s)" % (str(mediaStorage))
63         quit()
64
65     # Check Media Storage SOP Class
66     if (fileMetaInformation.FindDataElement(gdcm.Tag(0x0002, 0x0002))):
67         sopClassUid = str(fileMetaInformation.GetDataElement(gdcm.Tag(0x0002, 0x0002)).GetValue())
68         # Check SOP UID
69         if (sopClassUid != "1.2.840.10008.1.3.10"):
70             # File is not a DICOMDIR
71             print "This file is not a DICOMDIR"
72         else:
73             # Not present
```

```

73         print "Media Storage SOP Class not present"
74         quit()
75
76     # Iterate through the DICOMDIR data set
77     iterator = dataSet.GetDES().begin()
78     while (not iterator.equal(dataSet.GetDES().end())):
79         dataElement = iterator.next()
80
81         # Check the element tag
82         if (dataElement.GetTag() == gdcm.Tag(0x004, 0x1220)):
83             # The 'Directory Record Sequence' element
84             sequence = dataElement.GetValueAsSQ()
85
86             # Loop through the sequence items
87             itemNr = 1
88             while (itemNr < sequence.GetNumberOfItems()):
89                 item = sequence.GetItem(itemNr)
90
91                 # Check the element tag
92                 if (item.FindDataElement(gdcm.Tag(0x0004, 0x1430))):
93                     # The 'Directory Record Type' element
94                     value = str(item.GetDataElement(gdcm.Tag(0x0004, 0x1430)).GetValue())
95
96                     # PATIENT
97                     while (value.strip() == "PATIENT"):
98                         print value.strip()
99                         # Print patient name
100                        if (item.FindDataElement(gdcm.Tag(0x0010, 0x0010))):
101                            value = str(item.GetDataElement(gdcm.Tag(0x0010, 0x0010)).GetValue())
102                            print value
103
104                        # Print patient ID
105                        if (item.FindDataElement(gdcm.Tag(0x0010, 0x0020))):
106                            value = str(item.GetDataElement(gdcm.Tag(0x0010, 0x0020)).GetValue())
107                            print value
108
109                        # Next
110                        itemNr = itemNr + 1
111                        item = sequence.GetItem(itemNr)
112                        if (item.FindDataElement(gdcm.Tag(0x0004, 0x1430))):
113                            value = str(item.GetDataElement(gdcm.Tag(0x0004, 0x1430)).GetValue())
114
115                        # STUDY
116                        while (value.strip() == "STUDY"):
117                            print value.strip()
118
119                            # Print study UID
120                            if (item.FindDataElement(gdcm.Tag(0x0020, 0x000d))):
121                                value = str(item.GetDataElement(gdcm.Tag(0x0020, 0x000d)).GetValue())
122                                print value
123
124                            # Print study date
125                            if (item.FindDataElement(gdcm.Tag(0x0008, 0x0020))):
126                                value = str(item.GetDataElement(gdcm.Tag(0x0008, 0x0020)).GetValue())
127                                print value
128
129                            # Print study description
130                            if (item.FindDataElement(gdcm.Tag(0x0008, 0x1030))):
131                                value = str(item.GetDataElement(gdcm.Tag(0x0008, 0x1030)).GetValue())
132                                print value
133
134                            # Next
135                            itemNr = itemNr + 1
136                            item = sequence.GetItem(itemNr)
137                            if (item.FindDataElement(gdcm.Tag(0x0004, 0x1430))):
138                                value = str(item.GetDataElement(gdcm.Tag(0x0004, 0x1430)).
139                                GetValue())
140
141                            # SERIES
142                            while (value.strip() == "SERIES"):
143                                print value.strip()
144
145                                # Print series UID
146                                if (item.FindDataElement(gdcm.Tag(0x0020, 0x000e))):
147                                    value = str(item.GetDataElement(gdcm.Tag(0x0020, 0x000e)).
148                                    GetValue())
149                                    print value
150

```

```

149             # Print series modality
150             if (item.FindDataElement(gdcm.Tag(0x0008, 0x0060))):
151                 value = str(item.GetDataElement(gdcm.Tag(0x0008, 0x0060))).
152             GetValue()
153             print "Modality"
154             print value
155             # Print series description
156             if (item.FindDataElement(gdcm.Tag(0x0008, 0x103e))):
157                 value = str(item.GetDataElement(gdcm.Tag(0x0008, 0x103e))).
158             GetValue()
159             print "Description"
160             print value
161             # Next
162             itemNr = itemNr + 1
163             item = sequence.GetItem(itemNr)
164             if (item.FindDataElement(gdcm.Tag(0x0004, 0x1430))):
165                 value = str(item.GetDataElement(gdcm.Tag(0x0004, 0x1430))).
166             GetValue()
167             # IMAGE
168             while (value.strip() == "IMAGE"):
169                 print value.strip()
170             # Print image UID
171             if (item.FindDataElement(gdcm.Tag(0x0004, 0x1511))):
172                 value = str(item.GetDataElement(gdcm.Tag(0x0004, 0x1511))).
173             GetValue()
174             print value
175             # Next
176             if (itemNr < sequence.GetNumberOfItems()):
177                 itemNr = itemNr + 1
178             else:
179                 break
180             item = sequence.GetItem(itemNr)
181             if (item.FindDataElement(gdcm.Tag(0x0004, 0x1430))):
182                 value = str(item.GetDataElement(
183                     gdcm.Tag(0x0004, 0x1430)).GetValue())
184             # Next
185             itemNr = itemNr + 1

```

## 29.123 ReadAndPrintAttributes.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
/*
 * This small example will show how one can read and print
 * a DICOM Attribute using different technique (by tag or by name)
 */

#include "gdcmReader.h"
#include "gdcmGlobal.h"
#include "gdcmDicts.h"
#include "gdcmDict.h"
#include "gdcmAttribute.h"
#include "gdcmStringFilter.h"

#include <iostream>

int main(int argc, char *argv[])
{
    if( argc < 2 )

```

```

    {
        std::cerr << argv[0] << " input.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];

    // Instantiate the reader:
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Could not read: " << filename << std::endl;
        return 1;
    }

    // The output of gdcm::Reader is a gdcm::File
    gdcm::File &file = reader.GetFile();

    // the dataset is the the set of element we are interested in:
    gdcm::DataSet &ds = file.GetDataSet();

    const gdcm::Global& g = gdcm::Global::GetInstance();
    const gdcm::Dicts &dicts = g.GetDicts();
    const gdcm::Dict &pubdict = dicts.GetPublicDict();

    using namespace gdcm;

    // In this example we will show why using name to lookup attribute can be
    // dangerous.
    Tag tPatientName(0x00,0x00);
    //const DictEntry &de1 =
    pubdict.GetDictEntryByName("Patient Name", tPatientName);

    std::cout << "Found: " << tPatientName << std::endl;

    // Indeed the attribute could not be found. Since DICOM 2003, Patient Name
    // has become Patient's Name.

    Tag tPatientsName;
    //const DictEntry &de2 =
    pubdict.GetDictEntryByName("Patient's Name", tPatientsName);

    std::cout << "Found: " << tPatientsName << std::endl;

    // Let's try to read an arbitrary DICOM Attribute:
    Tag tDoseGridScaling;
    //const DictEntry &de3 =
    pubdict.GetDictEntryByName("Dose Grid Scaling", tDoseGridScaling);

    std::cout << "Found: " << tDoseGridScaling << std::endl;

    if( ds.FindDataElement( tDoseGridScaling ) )
    {
        gdcm::StringFilter sf;
        sf.SetFile(file);
        std::cout << "Attribute Value as String: " << sf.ToString( tDoseGridScaling ) << std::endl;

        // Let's check the name again:
        std::pair<std::string, std::string> pss
            = sf.ToStringPair( tDoseGridScaling );
        std::cout << "Attribute Name Checked: " << pss.first << std::endl;
        std::cout << "Attribute Value (string): " << pss.second << std::endl;

        //const DataElement &dgs = ds.GetDataElement( tDoseGridScaling );

        // Let's assume for a moment we knew the tag number:
        Attribute<0x3004,0x000e> at;
        assert( at.GetTag() == tDoseGridScaling );
        at.SetFromDataSet( ds );
        // For the sake of long term maintenance, we will not write
        // that this particular attribute is stored as a double. What if
        // a user made a mistake. It is much safer to rely on GDCM internal
        // mechanism to deduce the VR::DS type (represented as a iieee double)
        Attribute<0x3004,0x000e>::ArrayType v = at.
            GetValue();
        std::cout << "DoseGridScaling=" << v << std::endl;
    }

    return 0;
}

```

## 29.124 ReadExplicitLengthSQIVR.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmlReader.h"
#include "gdcmlImplicitDataElement.h"
#include "gdcmlDataSet.h"
#include "gdcmlPrivateTag.h"
#include "gdcmlPrivateTag.h"
#include "gdcmlByteValue.h"
#include "gdcmlSequenceOfItems.h"

using namespace gdcml;

int main(int argc, char *argv[])
{
    if ( argc < 2 ) return 1;
    const char *filename = argv[1];
    gdcml::Reader r;
    r.SetFileName( filename );
    r.Read();

    //gdcml::PrivateTag pt(0xel,0x42,"ELSCINT1");
    //gdcml::Tag pt(0x88,0x200);
    gdcml::Tag pt(0x8,0x1140);
    DataSet &ds = r.GetFile().GetDataSet();
    const DataElement &de = ds.GetDataElement( pt );

    std::cout << de << std::endl;
    const ByteValue *bv = de.GetByteValue();
    SmartPointer<SequenceOfItems> sqi = new
        SequenceOfItems;
    sqi->SetLength( bv->GetLength() );
    std::stringstream ss;
    ss.str( std::string( bv->GetPointer(), bv->GetLength() ) );
    sqi->Read<ImplicitDataElement,SwapperNoOp>( ss );

    std::cout << *sqi << std::endl;

    return 0;
}

```

## 29.125 ReadFiles.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
import gdcml.*;
import java.io.File;

public class ReadFiles
{
    static int i = 0;
    public static void process(String path)

```

```

{
    //String path = file.getPath();
    assert PosixEmulation.FileExists(path) : "Problem converting to 8bits";

    System.out.println("Reading: " + path );
    System.out.println("File: " + i++);
    Reader r = new Reader();
    try
    {
        r.SetFileName( path );
        TagSetType skip = new TagSetType();
        skip.insert( new Tag(0x7fe0,0x10) );
        boolean b = r.ReadUpToTag( new Tag(0x88,0x200), skip );
        //System.out.println("DS:\n" + r.GetFile().GetDataSet().toString() );
    }
    finally
    {
        r.delete(); // will properly call C++ destructor and close file descriptor
    }
}

// Process only files under dir
public static void visitAllFiles(File dir)
{
    if (dir.isDirectory())
    {
        String[] children = dir.list();
        for (int i=0; i<children.length; i++)
        {
            visitAllFiles(new File(dir, children[i]));
        }
    }
    else
    {
        process(dir.getPath());
    }
}

public static void waiting (int n)
{
    long t0, t1;
    t0 = System.currentTimeMillis();
    do
    {
        t1 = System.currentTimeMillis();
    }
    while ((t1 - t0) < (n * 1000));
}

public static void main(String[] args) throws Exception
{
    String directory = args[0];

    Directory gdir = new Directory();
    long n = gdir.Load( directory, true );
    System.out.println( gdir.toString() );
    FilenamesType files = gdir.GetFilenames();
    for( long i = 0; i < n; ++i )
    {
        String path = files.get( (int)i );
        process( path );
    }

    System.out.println( "Java API" );

    //waiting( 10 );
    for( int i = 0; i < 2; ++i )
    {
        File dir = new File(directory);
        visitAllFiles(dir);
    }
}
}

```

## 29.126 ReadGEMSSDO.cxx

```

/*=====

```



```

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmDataElement.h"
#include "gdcmPrivateTag.h"

#include <iostream>
#include <string>

using namespace gdcm;

struct SDOElement
{
    typedef std::vector<std::string>::size_type SizeType;
    const char *GetData(SizeType index) const {
        return Data[index].c_str();
    }
    SizeType GetNumberOfData() const {
        return Data.size();
    }
    void SetData(SizeType index, const char *data) {
        Data[index] = data;
    }
    const char *GetDataFormat() const {
        return DataFormat.c_str();
    }
    void SetDataFormat(const char *dataformat, SizeType num) {
        DataFormat = dataformat;
        Data.resize( num );
    }
    void Print( std::ostream &os ) const {
        os << DataFormat << ":" << std::endl;
        std::vector<std::string>::const_iterator it = Data.begin();
        size_t s = 0;
        for( ; it != Data.end(); ++it )
        {
            os << "    (" << s++ << ") " << *it << std::endl;
        }
    }
private:
    std::string DataFormat;
    std::vector<std::string> Data;
};

class SDOHeader
{
public:
    typedef std::vector<SDOElement> SDOElements;
    typedef SDOElements::size_type SizeType;
    SizeType GetNumberOfSDOElements() const {
        return InternalSDODataset.size();
    }
    void AddSDOElement(SDOElement const &sdoelement) {
        InternalSDODataset.push_back( sdoelement );
    }
    const SDOElement &GetSDOElement(SizeType index) const {
        return InternalSDODataset[index];
    }
    const SDOElement &GetSDOElementByName(const char *) const {
        return InternalSDODataset[0];
    }
    void LoadFromAttributes(std::string const &s1, std::string const &s2)
    {
        std::string tok;
        std::string tok2;
        std::stringstream strstr(s1);
        std::stringstream strstr2(s2);

        SDOElement element;
        // Do format
        size_t count = 0;

```

```

while ( std::getline ( strstr2, tok, '\\') )
{
    //std::cout << tok << " ";
    std::getline ( strstr2, tok2, '\\');
    //std::cout << tok2 << std::endl;
    count += atoi( tok2.c_str() );
    element.SetDataFormat( tok.c_str(), atoi( tok2.c_str() ) );
    for( size_t t = 0; t < element.GetNumberOfData(); ++t )
    {
        std::getline ( strstr, tok, '\\');
        element.SetData(t, tok.c_str() );
    }
    AddSDOElement( element );
}
//while ( std::getline ( strstr, tok, '^' ) )
// while ( std::getline ( strstr, tok, '\\') )
// {
//     std::cout << tok << std::endl;
//     count++;
// }
// std::cout << "Count: " << count << std::endl;
// count = 0;

// std::cout << "Count: " << count << std::endl;
}

void Print( std::ostream &os ) const {
    SDOElements::const_iterator it = InternalSDODataset.begin();
    for( ; it != InternalSDODataset.end(); ++it )
    {
        it->Print ( os );
    }
}

private:
    SDOElements InternalSDODataset;
};

bool sdo_decode( DataElement const &stringdata, DataElement const &stringdataformat )
{
    const char *sd = stringdata.GetByteValue()->GetPointer();
    const size_t len_sd = stringdata.GetByteValue()->GetLength();

    std::string s1 = std::string( sd, len_sd );

    const char *sdf = stringdataformat.GetByteValue()->GetPointer();
    const size_t len_sdf = stringdataformat.GetByteValue()->GetLength();

    std::string s2 = std::string( sdf, len_sdf );

    // std::cout << s1 << std::endl;
    // std::cout << s2 << std::endl;

    SDOHeader header;
    header.LoadFromAttributes( s1, s2 );

    header.Print( std::cout );

    return true;
}

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " input.dcm" << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        return 1;
    }

    File &file = reader.GetFile();
    DataSet &ds = file.GetDataSet();

    // StringData (0033,xx1F) 3 "GEMS_GENIE_1" List of SDO parameters stored as
    // list of strings
    const PrivateTag tstringdata(0x33,0x1f,"GEMS_GENIE_1");

```

```

// StringDataFormat (0033,xx23) 3 "GEMS_GENIE_1" Format of string parameters;
// contains information about name and number of strings in list
const PrivateTag tstringdataformat(0x33,0x23,"GEMS_GENIE_1");

if( !ds.FindDataElement( tstringdata ) ) return 1;
const DataElement& stringdata = ds.GetDataElement( tstringdata );
if( !ds.FindDataElement( tstringdataformat ) ) return 1;
const DataElement& stringdataformat = ds.GetDataElement( tstringdataformat );

sdo_decode( stringdata, stringdataformat );

return 0;
}

```

## 29.127 ReadMultiTimesException.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
// The intention of this sample program is to provoke bad_alloc exceptions in gdcm code

#include "gdcmImageReader.h"

int main(int argc, char* argv[])
{
    // We pre-allocate some memory (about 1Gb) to help the issue to show up earlier
    char *dummyBuffer = new char[1024*1024*1100]; (void)dummyBuffer;
    // Check the number of parameters given
    if (argc < 3)
    {
        std::cerr << "Usage: " << argv[0] << " Filename numberOfTries" << std::endl;
        return 1;
    }

    std::cout << "We are going to read the file: " << argv[1] << " " << argv[2] << " times" << std::endl;
    // We hold the pointers in an array to avoid the memory to be released
    // We read the input file n-times
    for (int i = 0; i < atoi(argv[2]); ++i)
    {
        gdcm::ImageReader reader;
        std::cout << "Reading try: " << i << std::endl;
        // Read files
        reader.SetFileName(argv[1]);
        try
        {
            reader.Read();
            gdcm::Image & img = reader.GetImage();
            unsigned long len = img.GetBufferLength();
            char *buffer = new char[ len ];
            img.GetBuffer( buffer ); // do NOT de-allocate buffer !
        }
        catch (std::bad_alloc)
        {
            std::cerr << "BAD ALLOC Exception caught!" << std::endl;
        }
        catch (...)
        {
            std::cerr << "Exception caught!" << std::endl;
        }
    }

    return 0;
}

```

## 29.128 ReadSeriesIntoVTK.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
// We are required to call the package 'vtk' eventhough I (MM) would have preferred
// an import statement along the line of:
// import vtkgdcm.*;
import vtk.*;

/*
 * Usage:
 * export LD_LIBRARY_PATH=/usr/lib/jvm/java-6-openjdk/jre/lib/amd64/xawt:.
 * java -classpath `pwd`/vtkgdcm.jar:/usr/share/java/vtk.jar:. ReadSeriesIntoVTK
 */
public class ReadSeriesIntoVTK
{
    static {
        System.loadLibrary("vtkCommonJava");
        System.loadLibrary("vtkFilteringJava");
        System.loadLibrary("vtkIOJava");
        System.loadLibrary("vtkImagingJava");
        System.loadLibrary("vtkGraphicsJava");
        System.loadLibrary("vtkgdcmJava");
        try {
            System.loadLibrary("vtkRenderingJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkHybrid, skipping...");
        }
        try {
            System.loadLibrary("vtkHybridJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkHybrid, skipping...");
        }
        try {
            System.loadLibrary("vtkVolumeRenderingJava");
        } catch (Throwable e) {
            System.out.println("cannot load vtkVolumeRendering, skipping...");
        }
    }

    public static void main(String[] args)
    {
        vtkFileOutputWindow outWin = new vtkFileOutputWindow();
        outWin.SetInstance(outWin);
        outWin.SetFileName("MVSvtkViewer.log");

        // See: http://review.source.kitware.com/#change,888
        // vtkWrapJava does not handle static keyword
        // String directory = vtkGDCMTesting.GetGDCMDataRoot();
        vtkGDCMTesting t = new vtkGDCMTesting();
        String directory = t.GetGDCMDataRoot();
        String file0 = directory + "/SIEMENS_MAGNETOM-12-MONO2-FileSeq0.dcm";
        String file1 = directory + "/SIEMENS_MAGNETOM-12-MONO2-FileSeq1.dcm";
        String file2 = directory + "/SIEMENS_MAGNETOM-12-MONO2-FileSeq2.dcm";
        String file3 = directory + "/SIEMENS_MAGNETOM-12-MONO2-FileSeq3.dcm";

        vtkStringArray s = new vtkStringArray();
        System.out.println("adding : " + file0 );
        s.InsertNextValue( file0 );
        s.InsertNextValue( file1 );
        s.InsertNextValue( file2 );
        s.InsertNextValue( file3 );

        vtkGDCMImageReader reader = new vtkGDCMImageReader();
        reader.SetFileNames( s );
        reader.Update();

        System.out.println("Success reading: " + file0 );
    }
}

```

```

    vtkMetaImageWriter writer = new vtkMetaImageWriter();
    writer.DebugOn();
    writer.SetCompression( false );
    writer.SetInput( reader.GetOutput() );
    writer.SetFileName( "ReadSeriesIntoVTK.mhd" );
    writer.Write();

    System.out.println("Success writing: " + writer.GetFileName() );
}
}

```

## 29.129 ReadUTF8QtDir.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * GDCM API expect a const char * as input for SetFileName
 * In order to use this API from Qt, here is a simple test that
 * shows how to do it in a portable manner:
 *
 * http://doc.qt.nokia.com/latest/qdir.html#navigation-and-directory-operations
 */

#include "gdcmlReader.h"
#include "gdcmlDirectory.h"

#include <QDir>
#include <QString>
#include <QCoreApplication>

#include <string>
#include <fstream>

#include <stdio.h> // fopen

static int TestBothFuncs(const char *info , const char *ba_str)
{
    int res = 0;
    FILE *f = fopen( ba_str, "r" );
    if( f )
    {
        std::cout << info << " fopen: " << ba_str << std::endl;
        fclose(f);
        ++res;
    }
    gdcml::Reader reader;
    std::ifstream is( ba_str, std::ios::binary );
    if( is.is_open() )
    {
        std::cout << info << " is_open: " << ba_str << std::endl;
        ++res;
    }
    reader.SetStream( is );
    if( reader.CanRead() == true )
    {
        std::cout << info << " SetStream/CanRead:" << ba_str << std::endl;
        ++res;
    }
    is.close();
    reader.SetFileName( ba_str );
    if( reader.CanRead() == true )
    {
        std::cout << info << " SetFileName/CanRead:" << ba_str << std::endl;
        ++res;
    }
    return 4 - res;
}

```

```

}

static int scanFolder(const char dirname[])
{
    int res = 0;
    gdcm::Directory dir;
    unsigned int nfiles = dir.Load( dirname, true );
    const gdcm::Directory::FileNamesType &filenames = dir.
        GetFileNames();

    for( unsigned int i = 0; i < nfiles; ++i )
    {
        const char *ba_str = filenames[i].c_str();
        res += TestBothFuncs("GDCM",ba_str);
    }
    return res;
}

static int scanFolderQt(QDir const &dir, QStringList& files)
{
    int res = 0;
    QFileInfoList children = dir.entryInfoList(QDir::AllEntries|QDir::NoDotAndDotDot);
    for ( int i=0; i<children.count(); i++ ) {
        QFileInfo file = children.at(i);
        if ( file.isDir() == true ) {
            res += scanFolderQt(QDir(file.absoluteFilePath()), files);
            continue;
        }
        // Convert back from the internal representation to 8bits
        // toLocal8Bit() returns by copy. Need to store explicitey the QByteArray
        QByteArray str = file.absoluteFilePath().toLocal8Bit();
        const char *ba_str1 = str.constData();
        res += TestBothFuncs("QString", ba_str1);
    }
    return res;
}

int main(int argc, char *argv[])
{
    // very important:
    QCoreApplication qCoreApp( argc , argv );
    if( argc < 2 )
    {
        std::cerr << argv[0] << " dir " << std::endl;
        return 1;
    }

    int res = 0;
    const char *dirname = argv[1];
    res += scanFolder( dirname );

    QDir dir( QString::fromLocal8Bit(dirname) );
    QStringList files;
    res += scanFolderQt( dir, files);

    if( res )
        std::cerr << "Problem with UTF-8" << std::endl;
    else
        std::cerr << "Success with UTF-8" << std::endl;

    return res;
}

```

## 29.130 RefCounting.cs

```

/*=====

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

```

```

=====*/
using Kitware.VTK;
using Kitware.VTK.GDCM;

/*
 * this is not so much an example but simply a test to make sure ctor / dctor work as expected
 * and call the ::New and ->Delete() of VTK style.
 */
public class RefCounting
{
    public static int Main(string[] args)
    {
        {
            vtkGDCMTesting testing1 = vtkGDCMTesting.New();
            vtkGDCMTesting testing2 = new vtkGDCMTesting(); // just in case people do
                not read STYLE documentation

            vtkGDCMImageReader reader1 = vtkGDCMImageReader.
                New();
            vtkGDCMImageReader reader2 = new vtkGDCMImageReader();

            vtkGDCMImageWriter writer1 = vtkGDCMImageWriter.
                New();
            vtkGDCMImageWriter writer2 = new vtkGDCMImageWriter();

            using (vtkGDCMTesting testing3 = new vtkGDCMTesting())
            {
                System.Console.Write( "GetReferenceCount: " + testing1.GetReferenceCount() + "\n");
                System.Console.Write( "GetReferenceCount: " + testing2.GetReferenceCount() + "\n");
                System.Console.Write( "GetReferenceCount: " + testing3.GetReferenceCount() + "\n");
            }

            using (vtkGDCMImageReader reader3 = new vtkGDCMImageReader())
            {
                System.Console.Write( "GetReferenceCount: " + reader3.GetReferenceCount() + "\n");
            }

            using (vtkGDCMImageWriter writer3 = vtkGDCMImageWriter.
                New())
            {
                System.Console.Write( "GetReferenceCount: " + writer3.GetReferenceCount() + "\n");
            }

            // C# destructor will call ->Delete on all C++ object as expected.
            return 0;
        }
    }
}

```

## 29.131 ReformatFile.cs

This is a C++ example on how to use FileDerivation

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcms.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Simple C# example
 *
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/ReformatFile.exe input.dcm output.dcm
 */
using System;
using gdcm;

public class ReformatFile

```

```

{
    public static int Main(string[] args)
    {
        gdcm.FileMetaInformation.
            SetSourceApplicationEntityTitle( "My Reformat App" );

        // http://www.oid-info.com/get/1.3.6.1.4.17434
        string THERALYS_ORG_ROOT = "1.3.6.1.4.17434";
        gdcm.UIDGenerator.SetRoot( THERALYS_ORG_ROOT );
        System.Console.WriteLine( "Root dir is now: " + gdcm.UIDGenerator.
            GetRoot() );

        string filename = args[0];
        string outfilename = args[1];

        Reader reader = new Reader();
        reader.SetFileName( filename );
        if( !reader.Read() )
        {
            System.Console.WriteLine( "Could not read: " + filename );
            return 1;
        }

        UIDGenerator uid = new UIDGenerator(); // helper for uid generation
        FileDerivation fd = new FileDerivation();
        // For the purpose of this exercise we will pretend that this image is referencing
        // two source image (we need to generate fake UID for that).
        string ReferencedSOPClassUID = "1.2.840.10008.5.1.4.1.1.7"; // Secondary Capture
        fd.AddReference( ReferencedSOPClassUID, uid.Generate() );
        fd.AddReference( ReferencedSOPClassUID, uid.Generate() );

        // Again for the purpose of the exercise we will pretend that the image is a
        // multiplanar reformat (MPR):
        // CID 7202 Source Image Purposes of Reference
        // { "DCM",121322,"Source image for image processing operation"},
        fd.SetPurposeOfReferenceCodeSequenceCodeValue( 121322 );
        // CID 7203 Image Derivation
        // { "DCM",113072,"Multiplanar reformatting" },
        fd.SetDerivationCodeSequenceCodeValue( 113072 );
        fd.SetFile( reader.GetFile() );
        // If all Code Value are ok the filter will execute properly
        if( !fd.Derive() )
        {
            return 1;
        }

        gdcm.FileMetaInformation fmi = reader.GetFile().GetHeader();
        // The following three lines make sure to regenerate any value:
        fmi.Remove( new gdcm.Tag(0x0002,0x0012) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0013) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0016) );

        Writer writer = new Writer();
        writer.SetFileName( outfilename );
        writer.SetFile( fd.GetFile() );
        if( !writer.Write() )
        {
            System.Console.WriteLine( "Could not write: " + outfilename );
            return 1;
        }

        return 0;
    }
}

```

## 29.132 RemovePrivateTags.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even

```



```

10 #         the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #         PURPOSE.  See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 python RemovePrivateTags.py input.dcm output.dcm
19 """
20
21 import sys
22 import gdcm
23
24
25 if __name__ == "__main__":
26
27     file1 = sys.argv[1]
28     file2 = sys.argv[2]
29
30     # Instanciate the reader.
31     r = gdcm.Reader()
32     r.SetFileName( file1 )
33     if not r.Read():
34         sys.exit(1)
35
36     # Remove private tags
37     ano = gdcm.Anonymizer()
38     ano.SetFile( r.GetFile() )
39     if not ano.RemovePrivateTags():
40         sys.exit(1)
41
42     # Write DICOM file
43     w = gdcm.Writer()
44     w.SetFile( ano.GetFile() )
45     #w.CheckFileMetaInformationOff() # Do not attempt to check meta header
46     w.SetFileName( file2 )
47     if not w.Write():
48         sys.exit(1)
49
50     # It is usually a good idea to exit the script with an error, as gdcm does not remove partial (incorrect)
51     # DICOM file
52     # (application level)

```

## 29.133 RescaleImage.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE.  See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/DecompressImage.exe gdcmData/012345.002.050.dcm rescaled.dcm
 */
using System;
using gdcm;

public class DecompressImage
{
    public static int Main(string[] args)
    {
        {
            string file1 = args[0];
            ImageReader reader = new ImageReader();
            reader.SetFileName( file1 );
            bool ret = reader.Read();
            if( !ret )

```

```

    {
        return 1;
    }

    Image image = reader.GetImage();
    PixelFormat pixeltype = image.GetPixelFormat();

    Rescaler r = new Rescaler();
    r.SetIntercept( 0 );
    r.SetSlope( 1.2 );
    r.SetPixelFormat( pixeltype );
    PixelFormat outputpt = new PixelFormat( r.ComputeInterceptSlopePixelFormat() );

    System.Console.WriteLine( "pixeltype" );
    System.Console.WriteLine( pixeltype.ToString() );
    System.Console.WriteLine( "outputpt" );
    System.Console.WriteLine( outputpt.ToString() );

    uint len = image.GetBufferLength();
    short[] input = new short[ len / 2 ]; // sizeof(short) == 2
    image.GetArray( input );

    double[] output = new double[ len / 2 ];
    r.Rescale( output, input, len );

    // First Pixel is:
    System.Console.WriteLine( "Input:" );
    System.Console.WriteLine( input[0] );

    System.Console.WriteLine( "Output:" );
    System.Console.WriteLine( output[0] );

    return 0;
}

```

## 29.134 reslicesphere.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
//
// Load a DICOM series.
// Position a sphere within the volume.
// Allow the user to change between Axial, Sagittal, Coronal, and
// Oblique view of the images and move through the slices.
// The display should show the resliced image and the cross section
// of the sphere intersecting that plane.
//

/*
from Scott Johnson /Scott Johnson neuwave com/
to VTK /vtkusers vtk.org/
date Tue, May 11, 2010 at 7:01 PM
*/
#include <sstream>
#include <string>

#include <vtkDICOMImageReader.h>
#include <vtkStringArray.h>
#include <vtkDirectory.h>
#include <vtkImageThreshold.h>
#include <vtkImageShiftScale.h>
#include <vtkImageReslice.h>
#include <vtkRenderWindowInteractor.h>
#include <vtkImageViewer2.h>

```

```

#include <vtkSphereSource.h>
#include <vtkPolyDataMapper.h>
#include <vtkPlane.h>
#include <vtkCutter.h>
#include <vtkActor.h>
#include <vtkCommand.h>
#include <vtkSmartPointer.h>
#include <vtkMatrix4x4.h>
#include <vtkInteractorObserver.h>
#include <vtkProperty.h>
#include <vtkRenderer.h>
#include <vtkImageData.h>
#include <vtkImageActor.h>
#include "vtkTransformPolyDataFilter.h"
#include <vtkCamera.h>
#include <vtkMath.h>
#include <vtkTransform.h>
#include <vtkTextActor.h>
#include <vtkActor2D.h>
#include <vtkPolyDataMapper2D.h>
#include <vtkProperty2D.h>
#include <vtkGDCMImageReader.h>
#include <vtkImageChangeInformation.h>

#include "gdcmDirectory.h"
#include "gdcmTesting.h"
#include "gdcmIPPSorter.h"

// Change to match the path to find Raw_0.vti or provide
// the parameter when starting ResliceSphere.

const double sphereCenter[3]={74, 219, 70};

// Angles (0, 0, 0)
const double AxialMatrix[] = { 1.0, 0.0, 0.0, 0.0,
                               0.0, 1.0, 0.0, 0.0,
                               0.0, 0.0, 1.0, 0.0,
                               0.0, 0.0, 0.0, 1.0 };

// Angles (0, 90, 0)
const double SagittalMatrix[] = { 0.0, 0.0, 1.0, 0.0,
                                   0.0, 1.0, 0.0, 0.0,
                                   -1.0, 0.0, 0.0, 0.0,
                                   0.0, 0.0, 0.0, 1.0 };

// Angles (-90, 0, 0)
const double CoronalMatrix[] = { 1.0, 0.0, 0.0, 0.0,
                                  0.0, 0.0, 1.0, 0.0,
                                  0.0, -1.0, 0.0, 0.0,
                                  0.0, 0.0, 0.0, 1.0 };

// Angles (0, 90, 31)
const double ObliqueMatrix[] = { 0.0, -0.515038, 0.857167, 0.0,
                                   0.0, 0.857167, 0.515038, 0.0,
                                   -1.0, 0.0, 0.0, 0.0,
                                   0.0, 0.0, 0.0, 1.0 };

class ResliceRender;

// Class to handle key press events.
class KeyCallback : public vtkCommand
{
public:
    static KeyCallback* New()
    {
        return new KeyCallback();
    }

    void Execute(vtkObject* caller, unsigned long eventId, void *calldata);
    void SetCallbackData(ResliceRender* reslice);

protected:
    ResliceRender* _reslice;
};

class ResliceRender
{
public:
    typedef enum _ORIENTATION
    {
        AXIAL = 0,
        SAGITTAL = 1,
        CORONAL = 2,
        OBLIQUE = 3
    }

```

```

    } ORIENTATION;

ResliceRender()
{
    _orientation=AXIAL;
}

~ResliceRender()
{
    _transform->Delete();
    _reader->Delete();
    _reslice->Delete();
    _interactor->Delete();
    _imageViewer->Delete();

    _sphere->Delete();
    _sphereMapper->Delete();
    _sphereActor->Delete();

    _plane->Delete();
    _cutter->Delete();
    _polyTransform->Delete();
    _ROIMapper->Delete();
    _ROIActor->Delete();

    _annotation->Delete();
}

void CreatePipeline(const char* fileName)
{
    vtkProperty2D* props;

    //_reader=vtkXMLImageDataReader::New();
    //_reader->SetFileName(fileName);
    //_reader->Update();

    //_reader=qzDICOMImageReader::New();
    _reader=vtkGDCMImageReader::New();

    //vtkDirectory *d = vtkDirectory::New();
    //d->Open(fileName);
    //d->Print( std::cout );
    gdcmm::Directory d;
    d.Load(fileName);
    gdcmm::Directory::FileNamesType const &files = d.
    GetFileNames();

    gdcmm::IPPSorter s;
    s.SetComputeZSpacing( true );
    s.SetZSpacingTolerance( 1e-3 );
    bool b = s.Sort( files );
    if( !b )
    {
        std::cerr << "Failed to sort:" << fileName << std::endl;
        //return ;
    }
    //std::cout << "Sorting succeeded:" << std::endl;
    //s.Print( std::cout );

    //std::cout << "Found z-spacing:" << std::endl;
    //std::cout << s.GetZSpacing() << std::endl;
    double ippzspacing = s.GetZSpacing();

    const std::vector<std::string> & sorted = s.GetFileNames();
    vtkStringArray *vtkfiles = vtkStringArray::New();
    std::vector< std::string >::const_iterator it = sorted.begin();
    for( ; it != sorted.end(); ++it)
    {
        const std::string &f = *it;
        vtkfiles->InsertNextValue( f.c_str() );
    }

    //_reader->SetDirectoryName(fileName);
    //_reader->SetFileNames( d->GetFiles() );
    _reader->SetFileNames( vtkfiles );
    _reader->Update();

    const vtkFloatingPointType *spacing = _reader->GetOutput()->GetSpacing();

    vtkImageChangeInformation *v16 = vtkImageChangeInformation::New();
    #if (VTK_MAJOR_VERSION >= 6)

```

```

v16->SetInputConnection( _reader->GetOutputPort() );
#else
v16->SetInput( _reader->GetOutput() );
#endif
v16->SetOutputSpacing( spacing[0], spacing[1], ippzspacing );
v16->Update();

_threshold=vtkImageThreshold::New();
_threshold->ThresholdByUpper(-3024.0);
_threshold->ReplaceOutOn();
_threshold->SetOutValue(0.0);
_threshold->SetInputConnection(v16->GetOutputPort());

_shift=vtkImageShiftScale::New();
_shift->SetShift(0);
_shift->SetScale(1);
_shift->SetInputConnection(_threshold->GetOutputPort());

// Initialize the reslice with an axial orientation.
vtkSmartPointer<vtkMatrix4x4> matrix =
    vtkSmartPointer<vtkMatrix4x4>::New();
matrix->Identity();

_transform = vtkTransform::New();
_transform->SetMatrix(matrix);

_reslice = vtkImageReslice::New();
_reslice->SetOutputDimensionality(3);

// PROBLEM:
// The original intent was to connect the same transform
// to the vtkImageReslice and vtkTransformPolyDataFilter,
// but the resulting reslices appear different using the
// vtkTransform as opposed to explicitly setting the
// reslice axes via SetResliceAxes. Also, if the vtkTransform
// is connected and orientated other than axial, the extents
// don't seem to update resulting in VTK believing the slice
// is out of range.

_reslice->SetResliceTransform(_transform);
_reslice->SetResliceAxes(matrix);
_reslice->SetInputConnection(_reader->GetOutputPort());
_reslice->SetInputConnection(_shift->GetOutputPort());

// Create the sphere target shape.
_sphere=vtkSphereSource::New();
_sphere->SetRadius(7.0);
_sphere->SetThetaResolution(16);
_sphere->SetPhiResolution(16);
_sphere->SetCenter(sphereCenter[0], sphereCenter[1], sphereCenter[2]);

_sphereMapper=vtkPolyDataMapper::New();
_sphereMapper->SetInputConnection(_sphere->GetOutputPort());

_sphereActor=vtkActor::New();
_sphereActor->SetMapper(_sphereMapper);
_sphereActor->PickableOff();
_sphereActor->GetProperty()->SetColor(1.0, 0.0, 0.0);
_sphereActor->GetProperty()->SetEdgeColor(1.0, 0.0, 0.0);
_sphereActor->GetProperty()->SetDiffuseColor(1.0, 0.0, 0.0);
_sphereActor->SetVisibility(true);

// Create the cutting pipeline.
// This plane will be positioned in the original image coordinate system.
_plane = vtkPlane::New();
_plane->SetNormal(0.0, 0.0, 1.0);

_cutter = vtkCutter::New();
_cutter->SetInputConnection(_sphere->GetOutputPort());
_cutter->SetCutFunction(_plane);
_cutter->GenerateCutScalarsOn();
_cutter->SetValue(0, 0.5);

// The transform attached to _polyTransform should move the cut
// ROI into the resliced coordinate system, which should be the
// same as the coordinate system of the resliced images.
// PROBLEM: It doesn't.
_polyTransform = vtkTransformPolyDataFilter::New();
_polyTransform->SetTransform(_transform);
_polyTransform->SetInputConnection(_cutter->GetOutputPort());

```

```

        _ROIMapper = vtkPolyDataMapper2D::New();
        _ROIMapper->SetInputConnection(_polyTransform->GetOutputPort());

    vtkCoordinate* coordinate = vtkCoordinate::New();
    coordinate->SetCoordinateSystemToWorld();
    _ROIMapper->SetTransformCoordinate(coordinate);

    _ROIActor = vtkActor2D::New();
    _ROIActor->SetMapper(_ROIMapper);

    // Make sure the cut can be seen, especially the edges.
    props=_ROIActor->GetProperty();
    props->SetLineWidth(2);
    props->SetOpacity(1.0);
    //
    // props->EdgeVisibilityOn();
    // props->SetDiffuse(0.8);
    // props->SetSpecular(0.3);
    // props->SetSpecularPower(20);
    // props->SetRepresentationToSurface();
    // props->SetDiffuseColor(1.0, 0.0, 0.0);
    // props->SetEdgeColor(1.0, 0.0, 0.0);
    props->SetColor(1.0, 0.0, 0.0);

    _interactor = vtkRenderWindowInteractor::New();

    // Create the image viewer and add the actor with the cut ROI.
    _imageView = vtkImageViewer2::New();
    _imageView->SetupInteractor(_interactor);
    _imageView->SetSize(400, 400);
    _imageView->SetColorWindow(1024);
    _imageView->SetColorLevel(800);
    _imageView->SetInputConnection(_reslice->GetOutputPort());
    _imageView->GetImageActor()->SetOpacity(0.5);

    _annotation = vtkTextActor::New();
    _annotation->SetTextScaleModeToViewport();
    _imageView->GetRenderer()->AddActor(_annotation);

    // Add the cut shape actor to the renderer.
    _imageView->GetRenderer()->AddActor(_ROIActor);

    // Set up the key handler.
    vtkSmartPointer<KeyCallback> callback = vtkSmartPointer<KeyCallback>::New();
    callback->SetCallbackData(this);
    _interactor->AddObserver(vtkCommand::KeyPressEvent, callback);

    _interactor->Initialize();
}

void Start()
{
    _interactor->Start();
}

void ResetOrientation()
{
    vtkSmartPointer<vtkMatrix4x4> matrix =
        vtkSmartPointer<vtkMatrix4x4>::New();
    matrix->Identity();

    SetOrientation(matrix);
}

// Make sure the orientation of the vtkImageReslice and
// vtkTransform are in sync.
void SetOrientation(vtkMatrix4x4* matrix)
{
    _reslice->SetResliceAxes(matrix);
    _reslice->Update();

    vtkMatrix4x4* inverse = vtkMatrix4x4::New();
    vtkMatrix4x4::Invert(matrix, inverse);

    _transform->SetMatrix(inverse);
    _transform->Update();
}

// Set the current slice of the current view.
void SetSlice(int slice)
{

```

```

        std::stringstream posString;

        double    center[3];
        double    spacing[3];
        double    origin[3];
        double    point[4];
        double    newPoint[4];

        vtkImageData* imageData;
        int newSlice;

        // Try to make sure the extents of the reslice are updated.
        // PROBLEM: It doesn't seem to work when changing the orientation.
        imageData=vtkImageData::SafeDownCast(_reslice->GetOutput());
    #if (VTK_MAJOR_VERSION >= 6)
        assert(0);
    #else
        imageData->UpdateInformation();
    #endif

    // Let vtkImageViewer2 handle the slice limits.
    _imageView->SetSlice(slice);
    newSlice=GetSlice();

    imageData->GetCenter(center);
    imageData->GetSpacing(spacing);
    imageData->GetOrigin(origin);

    // Compute the position of the center of the slice based on the
    // spacing of the slices. The resliced axis will always
    // be the "Z" axis.
    point[0]=center[0];
    point[1]=center[1];
    point[2]=(newSlice * spacing[2]) + origin[2];
    point[3]=1.0;

    // Convert the coordinate from the reslice coordinate system to the
    // original image coordinate system.
    // PROBLEM: Logically this seems like it should have been multiplied
    // by the inverse to translate from the resliced coordinate system to
    // the original coordinate system. However, multiplying by the inverse
    // sticks the plane in the wrong place completely. Using the original
    // matrix at least gets the Z coordinate right.
    vtkMatrix4x4* matrix=_reslice->GetResliceAxes();
    vtkSmartPointer<vtkMatrix4x4> inverse =
        vtkSmartPointer<vtkMatrix4x4>::New();
    vtkMatrix4x4::Invert(matrix, inverse);

    matrix->MultiplyPoint(point, newPoint);
    _plane->SetOrigin(newPoint[0], newPoint[1], newPoint[2]);

    // Annotate the image.
    posString << "Position: (" << newPoint[0] << ", " << newPoint[1]
        << ", " << newPoint[2] << ") Slice: " << newSlice;
    _annotation->SetInput(posString.str());

    _imageView->Render();
}

int GetSlice()
{
    return _imageView->GetSlice();
}

// Set the orientation of the view.
void SetOrientation(ResliceRender::ORIENTATION orientation)
{
    vtkCamera* camera=_imageView->GetRenderer()->GetActiveCamera();

    double spacing[3];
    double origin[3];
    double point[4];
    double newPoint[4];
    double initialPosition;
    double xDirCosine[3];
    double yDirCosine[3];
    double zDirCosine[3];
    double normal[3];

    vtkImageData* imageData;

```

```

    vtkSmartPointer<vtkMatrix4x4> matrix =
        vtkSmartPointer<vtkMatrix4x4>::New();

    _orientation=orientation;

    // Reset ViewUp
    camera->SetViewUp(0.0, 1.0, 0.0);

    // Compute the cut plane position to the input coordinate system.
    imageData=vtkImageData::SafeDownCast(_reslice->GetInput());
#ifdef (VTK_MAJOR_VERSION >= 6)
    assert(0);
#else
    imageData->UpdateInformation();
#endif
    imageData->GetSpacing(spacing);
    imageData->GetOrigin(origin);

    point[0]=origin[0];
    point[1]=origin[1];
    point[2]=origin[2];
    point[3]=1.0;

    switch (_orientation)
    {
    case AXIAL:
        matrix->DeepCopy(AxialMatrix);
        initialPosition=sphereCenter[2];
        break;

    case CORONAL:
        matrix->DeepCopy(CoronalMatrix);
        initialPosition=sphereCenter[1];
        break;

    case SAGITTAL:
        matrix->DeepCopy(SagittalMatrix);
        initialPosition=sphereCenter[0];
        break;

    case OBLIQUE:
        matrix->DeepCopy(ObliqueMatrix);
        initialPosition=sphereCenter[2];
        break;
    }

    // Move the origin from the original image coordinate system to the
    // resliced image coordinate system.
    matrix->MultiplyPoint(point, newPoint);
    matrix->SetElement(0, 3, newPoint[0]);
    matrix->SetElement(1, 3, newPoint[1]);
    matrix->SetElement(2, 3, newPoint[2]);

    ResetOrientation();
    SetOrientation(matrix);

    // Compute the cutting plane normal and set it.
    // PROBLEM: If the transformation is connected rather than
    // using SetResliceAxes, the Direction Cosines do not reflect
    // the orientation of the vtkImageReslice.
    _reslice->GetResliceAxesDirectionCosines(xDirCosine, yDirCosine,
                                              zDirCosine);
    vtkMath::Cross(xDirCosine, yDirCosine, normal);
    _plane->SetNormal(normal);

    // Set the extents and spacing of the reslice to account for
    // all of the data.
    _reslice->SetOutputExtentToDefault();
    _reslice->SetOutputSpacing(spacing[0], spacing[0], spacing[0]);

    // Force the vtkImageViewer2 to update.
    // PROBLEM: The whole extent does not seem to be set in time
    // for the first render. This results in an error because the
    // slice is positioned outside the old bounds.
#ifdef (VTK_MAJOR_VERSION >= 6)
    _imageView->SetInputData(NULL);
#else
    _imageView->SetInput(NULL);
#endif
    _imageView->SetInputConnection(_reslice->GetOutputPort());

```



```

        _imageView->GetRenderer()->ResetCameraClippingRange();
        _imageView->GetRenderer()->ResetCamera();

        // Set the initial slice to be at the center of the sphere.
        // Divide by the spacing because this will be undone in SetSlice.
        SetSlice( (int)(initialPosition / spacing[0]));
    }

    vtkRenderWindowInteractor* GetInteractor()
    {
        return _interactor;
    }

protected:
    ORIENTATION                _orientation;

    //qzDICOMImageReader*      _reader;
    vtkGDCMImageReader*        _reader;
    vtkImageThreshold*         _threshold;
    vtkImageShiftScale*        _shift;
    vtkImageReslice*           _reslice;
    vtkRenderWindowInteractor* _interactor;
    vtkImageViewer2*           _imageView;

    vtkSphereSource*           _sphere;
    vtkPolyDataMapper*          _sphereMapper;
    vtkActor*                   _sphereActor;

    vtkPlane*                   _plane;
    vtkCutter*                  _cutter;
    vtkTransform*               _transform;
    vtkTransformPolyDataFilter* _polyTransform;
    vtkPolyDataMapper2D*         _ROIMapper;
    vtkActor2D*                 _ROIActor;

    vtkTextActor*               _annotation;
};

// Catch KeyPress events.
// Up Arrow - increases the slice
// Down Arrow - decreases the slice
// 'A' - sets the view to Axial
// 'S' - sets the view to Sagittal
// 'C' - sets the view to Coronal
// 'O' - set the view to Oblique

void KeyCallback::Execute(vtkObject* caller, unsigned long eventId, void *calldata)
{
    (void)caller;
    (void)eventId;
    (void)calldata;
    std::string sym=_reslice->GetInteractor()->GetKeySym();

    if (!sym.compare("Up"))
    {
        _reslice->SetSlice(_reslice->GetSlice() + 1);
    }
    else if (!sym.compare("Down"))
    {
        _reslice->SetSlice(_reslice->GetSlice() - 1);
    }
    else if ((!sym.compare("A")) || (!sym.compare("a")))
    {
        _reslice->SetOrientation(ResliceRender::AXIAL);
    }
    else if ((!sym.compare("C")) || (!sym.compare("c")))
    {
        _reslice->SetOrientation(ResliceRender::CORONAL);
    }
    else if ((!sym.compare("S")) || (!sym.compare("s")))
    {
        _reslice->SetOrientation(ResliceRender::SAGITTAL);
    }
    else if ((!sym.compare("O")) || (!sym.compare("o")))
    {
        _reslice->SetOrientation(ResliceRender::OBLIQUE);
    }
}

void KeyCallback::SetCallbackData(ResliceRender* reslice)

```

```

{
    _reslice=reslice;
}

// Usage: ResliceSphere [fileName]
int main(int argc, char *argv[])
{
    ResliceRender render;

    if (argc == 1)
    {
        const char *root = gdcm::Testing::GetDataExtraRoot();
        std::string dir3 = root;
        dir3 += "/gdcmSampleData/ForSeriesTesting/Dentist/images/";
        render.CreatePipeline(dir3.c_str());
    }
    else
    {
        render.CreatePipeline(argv[1]);
    }

    render.SetOrientation(ResliceRender::AXIAL);
    render.Start();

    return EXIT_SUCCESS;
}

```

## 29.135 ReWriteSCAsMR.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 GDCM 1.x would write out MR Image Storage as Secondary Capture Object while still setting Rescale
17 Slope/Intercept
18 and saving the Pixel Spacing in (0028,0030)
19 """
20 import gdcm
21 import sys,os
22
23 def CheckSecondaryCaptureObjectIsMRImageStorage(r):
24     ds = r.GetFile().GetDataSet()
25     # Check Source Image Sequence
26     if ds.FindDataElement( gdcm.Tag(0x0008,0x2112) ):
27         sis = ds.GetDataElement( gdcm.Tag(0x0008,0x2112) )
28         sqsis = sis.GetSequenceOfItems()
29         if sqsis.GetNumberOfItems():
30             item1 = sqsis.GetItem(1)
31             nestedds = item1.GetNestedDataSet()
32             if nestedds.FindDataElement( gdcm.Tag(0x0008,0x1150) ):
33                 ReferencedSOPClassUID = nestedds.GetDataElement( gdcm.Tag(0x0008,0x1150) )
34                 raw = ReferencedSOPClassUID.GetByteValue().GetPointer()
35                 uids = gdcm.UIDs()
36                 # what is the actual object we are looking at ?
37                 ms = gdcm.MediaStorage()
38                 ms.SetFromDataSet(ds)
39                 msuid = ms.GetString()
40                 uids.SetFromUID( msuid )
41                 msuidname = uids.GetName() # real Media Storage Name
42                 uids.SetFromUID( raw )
43                 sqmsuidname = uids.GetName() # Source Image Sequence Media Storage Name
44                 # If object is SC and Source derivation is MRImageStorage then we can assume 'Pixel Spacing' is
45                 correct
46                 if( sqmsuidname == 'MR Image Storage' and msuidname == 'Secondary Capture Image Storage' ):
47                     return True

```

```

47  # in all other case simply return the currentspacing:
48  return False
49
50  if __name__ == "__main__":
51      r = gdcm.ImageReader()
52      filename = sys.argv[1]
53      r.SetFileName( filename )
54      if not r.Read():
55          sys.exit(1)
56      f = r.GetFile()
57
58      if( CheckSecondaryCaptureObjectIsMRImageStorage(r) ):
59          # Special handling of the spacing:
60          # GDCM 1.2.0 would not rewrite correctly DICOM Object and would always set them as 'Secondary Capture
        Image Storage'
61          # while we would rather have 'MR Image Storage'
62          gdcm.ImageHelper.SetForcePixelSpacing( True )
63          mrspacing = gdcm.ImageHelper.GetSpacingValue( r.GetFile() )
64          # TODO: I cannot do simply the following:
65          #image.SetSpacing( mrspacing )
66          image.SetSpacing(0, mrspacing[0] )
67          image.SetSpacing(1, mrspacing[1] )
68          image.SetSpacing(2, mrspacing[2] )
69          gdcm.ImageHelper.SetForceRescaleInterceptSlope( True )
70          ris = gdcm.ImageHelper.GetRescaleInterceptSlopeValue(
            r.GetFile() )
71          image.SetIntercept( ris[0] )
72          image.SetSlope( ris[1] )
73
74      outfilename = sys.argv[2]
75      w = gdcm.ImageWriter()
76      w.SetFileName( outfilename )
77      w.SetFile( r.GetFile() )
78      w.SetImage( image )
79      if not w.Write():
80          sys.exit(1)
81
82      sys.exit(0)

```

## 29.136 rle2img.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
 * This example shows how to rewrite a ELSCINT1/PMSCT_RLE1 compressed
 * image so that it is readable by most 3rd party software (DICOM does
 * not specify this particular encoding).
 * This is required for the sake of interoperability with any standard
 * conforming DICOM system.
 *
 * Everything done in this code is for the sole purpose of writing interoperable
 * software under Sect. 1201 (f) Reverse Engineering exception of the DMCA.
 * If you believe anything in this code violates any law or any of your rights,
 * please contact us (gdcm-developers@lists.sourceforge.net) so that we can
 * find a solution.
 *
 * Everything you do with this code is at your own risk, since decompression
 * algorithm was not written from specification documents.
 *
 * Special thanks to:
 * Mauro Maiorca for bringing to our attention on this new ELSCINT1
 * compression algorithm : PMSCT_RLE1 (different from the 'LOSSLESS RICE')
 * See post at:
 * http://groups.google.com/group/comp.protocols.dicom/msg/f2b99bf706a7f8ca
 *
 * Thanks to Jesus Spinola, for more datasets,

```

```

* http://www.itk.org/pipermail/insight-users/2008-April/025571.html
*
* And last but not least, a very big thank to Ivo van Poorten, without
* whom we would still be looking at this compressed byte stream as if
* it was RLE compressed.
*/
#include "gdcmReader.h"
#include "gdcmPrivateTag.h"
#include "gdcmAttribute.h"
#include "gdcmImageWriter.h"

/* FIXME: Why is PhilipsLosslessRice.dcm a 512x512 image ... */
void delta_decode(const char *inbuffer, size_t length, std::vector<unsigned short> &output)
{
    // RLE pass
    std::vector<char> temp;
    for(size_t i = 0; i < length; ++i)
    {
        if( inbuffer[i] == (char)0xa5 )
        {
            //unsigned char repeat = (unsigned char)inbuffer[i+1] + 1;
            //assert( (unsigned char)inbuffer[i+1] != 255 );
            int repeat = (unsigned char)inbuffer[i+1] + 1;
            char value = inbuffer[i+2];
            while(repeat)
            {
                temp.push_back( value );
                --repeat;
            }
            i+=2;
        }
        else
        {
            temp.push_back( inbuffer[i] );
        }
    }

    // Delta encoding pass
    unsigned short delta = 0;
    for(size_t i = 0; i < temp.size(); ++i)
    {
        if( temp[i] == 0xa5 )
        {
            unsigned char v1 = (unsigned char)temp[i+1];
            unsigned char v2 = (unsigned char)temp[i+2];
            unsigned short value = (unsigned short)(v2 * 256 + v1);
            output.push_back( value );
            delta = value;
            i+=2;
        }
        else
        {
            unsigned short value = (unsigned short)(temp[i] + delta);
            output.push_back( value );
            delta = value;
        }
        //assert( output[output.size()-1] == ref[output.size()-1] );
    }

    if ( output.size() % 2 )
    {
        output.resize( output.size() - 1 );
    }
    std::cout << length << " -> " << output.size() * 2 << std::endl;
}

int main(int argc, char *argv [])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << "input.dcm [output.dcm]" << std::endl;
        std::cerr << "will default to 'outrle.dcm' unless output.dcm is specified."
        << std::endl;
        return 1;
    }
    const char *filename = argv[1];
    gdcm::Reader reader;
    reader.SetFileName( filename );
    if( !reader.Read() )
    {
        std::cerr << "Failed to read: " << filename << std::endl;
    }
}

```

```

    return 1;
}
const gdcm::DataSet& ds = reader.GetFile().GetDataSet();

// (07a1,1011) CS [PMSCT_RLE1] # 10,1 Tamar Compression Type
const gdcm::PrivateTag tcompressiontype(0x07a1,0x0011,"ELSCINT1");
if( !ds.FindDataElement( tcompressiontype ) ) return 1;
const gdcm::DataElement& compressiontype = ds.GetDataElement(
    tcompressiontype );
if ( compressiontype.IsEmpty() ) return 1;
const gdcm::ByteValue * bv = compressiontype.GetByteValue();
std::string comprle = "PMSCT_RLE1";
std::string comprgb = "PMSCT_RGB1";
bool isrle = false;
bool isrgb = false;
if( strcmp( bv->GetPointer(), comprle.c_str(), comprle.size() ) == 0 )
{
    isrle = true;
}
if( strcmp( bv->GetPointer(), comprgb.c_str(), comprgb.size() ) == 0 )
{
    isrgb = true;
    std::cerr << "See: pmsct_rgb1.cxx instead" << std::endl;
    return 1;
}
if( !isrgb && !isrle ) return 1;

const gdcm::PrivateTag tcompressedpixeldata(0x07a1,0x000a,"ELSCINT1");
if( !ds.FindDataElement( tcompressedpixeldata ) ) return 1;
const gdcm::DataElement& compressionpixeldata = ds.
    GetDataElement( tcompressedpixeldata );
if ( compressionpixeldata.IsEmpty() ) return 1;
const gdcm::ByteValue * bv2 = compressionpixeldata.GetByteValue();

gdcm::Attribute<0x0028,0x0010> at1;
at1.SetFromDataSet( ds );
gdcm::Attribute<0x0028,0x0011> at2;
at2.SetFromDataSet( ds );

gdcm::DataElement pixeldata( gdcm::Tag(0x7fe0,0x0010) );
pixeldata.SetVR( gdcm::VR::OW );
gdcm::VL bv2l = bv2->GetLength();
gdcm::VL at1l = at1.GetValue() * at2.GetValue() * 2; /* sizeof(unsigned short) ==
    2 */
// Handle special case that is not compressed:
if( bv2l == at1l )
{
    pixeldata.SetByteValue( bv2->GetPointer(), bv2->GetLength() );
}
else
{
    std::vector<unsigned short> buffer;
    delta_decode(bv2->GetPointer(), bv2->GetLength(), buffer);
    pixeldata.SetByteValue( (char*)&buffer[0], (uint32_t)(buffer.size() * sizeof( unsigned short )) );
}
// TODO we should check that decompress byte buffer match the expected size (row*col*...)

// Add the pixel data element
reader.GetFile().GetDataSet().Replace( pixeldata );

reader.GetFile().GetHeader().SetDataSetTransferSyntax(
    gdcm::TransferSyntax::ExplicitVRLittleEndian);
gdcm::Writer writer;
writer.SetFile( reader.GetFile() );

// Cleanup stuff:
// remove the compressed pixel data:
// FIXME: should I remove more private tags ? all of them ?
// oh well this is just an example
// use gdcm::Anonymizer::RemovePrivateTags if needed...
writer.GetFile().GetDataSet().Remove( compressionpixeldata.
    GetTag() );
std::string outfilename;
if (argc > 2)
    outfilename = argv[2];
else
    outfilename = "outrle.dcm";
writer.SetFileName( outfilename.c_str() );
if( !writer.Write() )
{
    std::cerr << "Failed to write" << std::endl;
}

```

```

    return 1;
}

std::cout << "success !" << std::endl;

return 0;
}

```

## 29.137 rtstructapp.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "vtkGDCMPolyDataReader.h"
#include "vtkGDCMPolyDataWriter.h"

#include "vtkPolyDataWriter.h"
#include "vtkPolyDataMapper.h"
#include "vtkPolyDataMapper2D.h"
#include "vtkActor2D.h"
#include "vtkRenderWindowInteractor.h"
#include "vtkMedicalImageProperties.h"
#include "vtkRenderWindow.h"
#include "vtkRenderer.h"
#include "vtkCamera.h"
#include "vtkProperty.h"
#include "vtkProperty2D.h"
#include "vtkAppendPolyData.h"
#include "vtkImageData.h"

/*
 * Small example to read in a RTSTRUCT and write it out (displays it too).
 */

// gdcmDataExtra/gdcmlNonImageData/exRT_Structure_Set_Storage.dcm
// gdcmDataExtra/gdcmlNonImageData/RTSTRUCT_1.3.6.1.4.1.22213.1.1396.2.dcm
// gdcmDataExtra/gdcmlNonImageData/RT/RTStruct.dcm

int main(int argc, char *argv[])
{
    if( argc < 3 )
    {
        std::cerr << argv[0] << " input.dcm output.dcm\n";
        return 1;
    }
    const char * filename = argv[1];
    const char * outfilename = argv[2];
    vtkGDCMPolyDataReader * reader =
        vtkGDCMPolyDataReader::New();
    reader->SetFileName( filename );
    reader->Update();

    //std::cout << reader->GetMedicalImageProperties()->GetStudyDate() << std::endl;

    vtkGDCMPolyDataWriter * writer =
        vtkGDCMPolyDataWriter::New();
    writer->SetNumberOfInputPorts( reader->GetNumberOfOutputPorts() );
    writer->SetFileName( outfilename );
    for(int num = 0; num < reader->GetNumberOfOutputPorts(); ++num )
    #if (VTK_MAJOR_VERSION >= 6)
        writer->SetInputConnection( num, reader->GetOutputPort(num) );
    #else
        writer->SetInput( num, reader->GetOutput(num) );
    #endif
    //doesn't look like the medical properties are actually written out
    writer->SetMedicalImageProperties( reader->GetMedicalImageProperties() );
    writer->SetRTStructSetProperties( reader->GetRTStructSetProperties() );
}

```

```

writer->Write();

// print reader output:
reader->Print( std::cout );
// print first output:
reader->GetOutput()->Print( std::cout );

vtkAppendPolyData *append = vtkAppendPolyData::New();

int n = reader->GetNumberOfOutputPorts();
for(int i = 0; i < n; ++i)
{
#ifdef VTK_MAJOR_VERSION >= 6
    append->AddInputConnection( reader->GetOutputPort(i) );
#else
    append->AddInput( reader->GetOutput(i) );
#endif
}

// Now we'll look at it.
vtkPolyDataMapper *cubeMapper = vtkPolyDataMapper::New();
#ifdef VTK_MAJOR_VERSION >= 6
    cubeMapper->SetInputConnection( append->GetOutputPort() );
#else
    cubeMapper->SetInput( append->GetOutput() );
#endif
cubeMapper->SetScalarRange(0,7);
vtkActor *cubeActor = vtkActor::New();
cubeActor->SetMapper(cubeMapper);
vtkProperty *property = cubeActor->GetProperty();
property->SetRepresentationToWireframe();

vtkRenderer *renderer = vtkRenderer::New();
vtkRenderWindow *renWin = vtkRenderWindow::New();
renWin->AddRenderer(renderer);

vtkRenderWindowInteractor *iren = vtkRenderWindowInteractor::New();
iren->SetRenderWindow(renWin);

renderer->AddActor(cubeActor);
renderer->ResetCamera();
renderer->SetBackground(1,1,1);

renWin->SetSize(300,300);

renWin->Render();
iren->Start();

reader->Delete();
append->Delete();
cubeMapper->Delete();
cubeActor->Delete();
renderer->Delete();
renWin->Delete();
iren->Delete();
writer->Delete();

return 0;
}

```

## 29.138 ScanDirectory.cs

This is a C# example on how to use Scanner

```

/*=====

Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====

```

```

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/ScanDirectory.exe /path/to/gdcmData/
 */
using System;
using gdcm;

public class ScanDirectory
{
    public static int Main(string[] args)
    {
        string directory = args[0];
        Tag t = new Tag(0x8,0x8);

        Directory d = new Directory();
        uint nfiles = d.Load( directory );
        if(nfiles == 0) return 1;
        //System.Console.WriteLine( "Files:\n" + d.toString() );

        //Scanner s = new Scanner();
        SmartPtrScan sscan = Scanner.New();
        Scanner s = sscan.__ref__();
        SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(s, "MySimple");
        s.AddTag( t );
        bool b = s.Scan( d.GetFileNames() );
        if(!b) return 1;

        System.Console.WriteLine( "Scan:\n" + s.toString() );

        System.Console.WriteLine( "success" );
        return 0;
    }
}

```

## 29.139 ScanDirectory.java

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

import gdcm.*;
import gdcm.Reader;
import gdcm.LookupTable;
import java.io.File;
import java.io.*;
import java.awt.image.*;
import javax.imageio.ImageIO;

public class ScanDirectory
{
    public static class MyWatcher extends SimpleSubjectWatcher
    {
        public MyWatcher(Subject s) { super(s,"Override String"); }
        protected void ShowProgress(Subject caller, Event evt)
        {
            ProgressEvent pe = ProgressEvent.Cast(evt);
            System.out.println( "This is my progress: " + pe.GetProgress() );
        }
    }

    public static byte[] GetAsByte(Bitmap input)
    {
        long len = input.GetBufferLength();
        byte[] buffer = new byte[ (int)len ];
    }
}

```



```

PhotometricInterpretation pi = input.GetPhotometricInterpretation();
if( pi.GetType() == PhotometricInterpretation.PIType.MONOCHROME1 )
{
    ImageChangePhotometricInterpretation icpi = new ImageChangePhotometricInterpretation();
    icpi.SetInput( input );
    icpi.SetPhotometricInterpretation(
        new PhotometricInterpretation(
            PhotometricInterpretation.PIType.MONOCHROME2 ) );
    if( icpi.Change() )
    {
        Bitmap output = icpi.GetOutput();
        output.GetArray( buffer );
    }
    return buffer;
}
else
{
    input.GetArray( buffer );
    return buffer;
}
}

public static short[] GetAsShort(Bitmap input)
{
    long len = input.GetBufferLength(); // length in bytes
    short[] buffer = new short[ (int)len / 2 ];
    PhotometricInterpretation pi = input.GetPhotometricInterpretation();
    if( pi.GetType() == PhotometricInterpretation.PIType.MONOCHROME1 )
    {
        ImageChangePhotometricInterpretation icpi = new ImageChangePhotometricInterpretation();
        icpi.SetInput( input );
        icpi.SetPhotometricInterpretation(
            new PhotometricInterpretation(
                PhotometricInterpretation.PIType.MONOCHROME2 ) );
        if( icpi.Change() )
        {
            Bitmap output = icpi.GetOutput();
            output.GetArray( buffer );
        }
        return buffer;
    }
    else
    {
        input.GetArray( buffer );
        return buffer;
    }
}

public static boolean WritePNG(Bitmap input, String outfilename )
{
    int imageType = BufferedImage.TYPE_CUSTOM;
    PixelFormat pf = input.GetPixelFormat();
    PhotometricInterpretation pi = input.GetPhotometricInterpretation();
    // We need to handle both public and private icon
    // It could well be that we are getting an RGB Icon or 16 bits Icon:
    ColorModel colorModel = null;
    if( pf.GetSamplesPerPixel() == 1 )
    {
        if( pi.GetType() == PhotometricInterpretation.PIType.MONOCHROME1
            || pi.GetType() == PhotometricInterpretation.PIType.MONOCHROME2 )
        {
            if( pf.GetScalarType() == PixelFormat.ScalarType.UINT8 )
            {
                imageType = BufferedImage.TYPE_BYTE_GRAY;
            }
            else if( pf.GetScalarType() == PixelFormat.ScalarType.UINT12 )
            {
                imageType = BufferedImage.TYPE_USHORT_GRAY;
            }
            else if( pf.GetScalarType() == PixelFormat.ScalarType.UINT16 )
            {
                imageType = BufferedImage.TYPE_USHORT_GRAY;
            }
        }
        else if( pi.GetType() == PhotometricInterpretation.PIType.PALETTE_COLOR )
        {
            LookupTable lut = input.GetLUT();
            long r1 = lut.GetLUTLength( LookupTable.LookupTableType.RED );
            byte[] rbuf = new byte[ (int)r1 ];
            long r12 = lut.GetLUT( LookupTable.LookupTableType.RED, rbuf );
            assert r1 == r12;
            long g1 = lut.GetLUTLength( LookupTable.LookupTableType.GREEN );
            byte[] gbuf = new byte[ (int)g1 ];

```

```

        long gl2 = lut.GetLUT( LookupTable.LookupTableType.GREEN, gbuf );
        assert gl == gl2;
        long bl = lut.GetLUTLength( LookupTable.LookupTableType.BLUE );
        byte[] bbuf = new byte[ (int)bl ];
        long bl2 = lut.GetLUT( LookupTable.LookupTableType.BLUE, bbuf );
        assert bl == bl2;
        colorModel = new IndexColorModel(8, (int)rl, rbuf, gbuf, bbuf);
        // For code below
        imageType = BufferedImage.TYPE_BYTE_GRAY;
    }
}
else if( pf.GetSamplesPerPixel() == 3 )
{
    if( pf.GetScalarType() == PixelFormat.ScalarType.UINT8 )
    {
        // FIXME should be TYPE_3BYTE_RGB
        imageType = BufferedImage.TYPE_3BYTE_BGR;
    }
}
//System.out.println( "pf: " + pf.toString() );
//System.out.println( "pi: " + pi.toString() );
long width = input.GetDimension(0);
long height = input.GetDimension(0);
BufferedImage bi;
if( pi.GetType() == PhotometricInterpretation.PIType.PALETTE_COLOR )
{
    bi = new BufferedImage(colorModel,
        colorModel.createCompatibleWritableRaster((int)width, (int)height),
        false, null);
}
else
{
    bi = new BufferedImage((int)width, (int)height, imageType);
}
WritableRaster wr = bi.getRaster();
//System.out.println( "imagetype: " + imageType );
if( imageType == BufferedImage.TYPE_BYTE_GRAY
    || imageType == BufferedImage.TYPE_3BYTE_BGR )
{
    byte[] buffer = GetAsByte( input );
    wr.setDataElements (0, 0, (int)width, (int)height, buffer);
}
else if( imageType == BufferedImage.TYPE_USHORT_GRAY )
{
    short[] buffer = GetAsShort( input );
    wr.setDataElements (0, 0, (int)width, (int)height, buffer);
}

File outputfile = new File( outfilename );
try {
    ImageIO.write(bi, "png", outputfile);
} catch (IOException e) {
    return false;
}
return true;
}

public static void main(String[] args) throws Exception
{
    String directory = args[0];

    Directory d = new Directory();
    long nfiles = d.Load( directory, true );
    if(nfiles == 0)
    {
        throw new Exception("No files found");
    }
    // System.out.println( "Files:\n" + d.toString() );
    FilenamesType fns = d.GetFilenames();

    //Scanner s = new Scanner();
    SmartPtrScan sscan = Scanner.New();
    Scanner s = sscan.__ref__();
    //SimpleSubjectWatcher watcher = new SimpleSubjectWatcher(s, "MySimple");
    MyWatcher watcher = new MyWatcher(s);
    Tag[] tagarray = {
        new Tag(0x0010, 0x0010), // PatientName
        new Tag(0x0010, 0x0020), // PatientID
        new Tag(0x0010, 0x0030), // PatientBirthDate
        new Tag(0x0010, 0x0040), // PatientSex
        new Tag(0x0010, 0x1010), // PatientAge
    }

```

```

        new Tag(0x0020, 0x000d), // StudyInstanceUID
        new Tag(0x0020, 0x0010), // StudyID
        new Tag(0x0008, 0x0020), // StudyDate
        new Tag(0x0008, 0x1030), // StudyDescription
        new Tag(0x0020, 0x000e), // SeriesInstanceUID
        new Tag(0x0020, 0x0011), // SeriesNumber
        new Tag(0x0008, 0x0021), // SeriesDate
        new Tag(0x0008, 0x103e), // SeriesDescription
        new Tag(0x0008, 0x0090), // ReferringPhysicianName
        new Tag(0x0008, 0x0060), // Modality
        new Tag(0x0054, 0x0400), // ImageID ?? Should be Instance number ??
        new Tag(0x0008, 0x0018), // SOPInstanceUID
        new Tag(0x0008, 0x0032), // AcquisitionTime
        new Tag(0x0008, 0x0033), // ContentTime
        new Tag(0x0020, 0x0013), // InstanceNumber
        new Tag(0x0020, 0x1041), // SliceLocation
        new Tag(0x0018, 0x0050), // SliceThickness ?? Eg. Enhanced MR Image Storage
        new Tag(0x0008, 0x0080), // InstitutionName
        new Tag(0x0028, 0x1050), // WindowCenter
        new Tag(0x0028, 0x1051), // WindowWidth
    };
    for( Tag t : tagarray ) {
        //System.out.println( "Tag: " + t.toString() );
        s.AddTag( t );
    }
    boolean b = s.Scan( fns );
    if(!b)
    {
        throw new Exception("Could not scan");
    }

    for( long idx = 0; idx < fns.size(); ++idx )
    {
        Reader r = new Reader();
        String fn = fns.get( (int)idx );
        String outfn = fn + ".png";
        r.SetFileName( fn );
        TagSetType tst = new TagSetType();
        tst.insert( new Tag(0x7fe0,0x10) );
        b = r.ReadUpToTag( new Tag(0x88,0x200), tst );
        UIntArrayType dims = ImageHelper.GetDimensionsValue( r.GetFile() );
        if( b )
        {
            IconImageFilter iif = new IconImageFilter();
            System.out.println( "Processing: " + fn );

            iif.SetFile( r.GetFile() );
            b = iif.Extract();
            if( b )
            {
                Bitmap icon = iif.GetIconImage(0);
                WritePNG(icon, outfn);
            }
            else
            {
                ImageReader ir = new ImageReader();
                ir.SetFileName( fn );
                if( ir.Read() )
                {
                    Image img = ir.GetImage();
                    StringFilter sf = new StringFilter();
                    sf.SetFile( r.GetFile() );
                    String strval = sf.ToString( new Tag(0x0028,0x0120) );
                    IconImageGenerator iig = new IconImageGenerator();
                    iig.SetPixmap( img );
                    iig.AutoPixelMinMax( true );
                    try {
                        double val = Double.parseDouble( strval );
                        iig.SetOutsideValuePixel( val );
                    }
                    catch ( NumberFormatException e ) {
                        //
                    }
                    iig.ConvertRGBToPaletteColor( false );
                    long idims[] = { 128, 128 };
                    iig.SetOutputDimensions( idims );
                    iig.Generate();
                    Bitmap icon = iig.GetIconImage();
                    WritePNG(icon, outfn);
                }
            }
        }
    }
}

```

```

    }

    System.out.println( "Scan:\n" + s.toString() );

    System.out.println( "success" );
}
}

```

## 29.140 ScanDirectory.py

```

1 #####
2 #
3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 import gdcm
16 import sys,os
17
18 class ProgressWatcher(gdcm.SimpleSubjectWatcher):
19     def ShowProgress(self, sender, event):
20         pe = gdcm.ProgressEvent.Cast(event)
21         print pe.GetProgress()
22     def EndFilter(self):
23         print "Yay ! I am done"
24
25 if __name__ == "__main__":
26     directory = sys.argv[1]
27
28     # Define the set of tags we are interested in
29     t1 = gdcm.Tag(0x8,0x8);
30     t2 = gdcm.Tag(0x10,0x10);
31
32     # Iterate over directory
33     d = gdcm.Directory();
34     nfiles = d.Load( directory );
35     if(nfiles == 0): sys.exit(1);
36     # System.Console.WriteLine( "Files:\n" + d.toString() );
37
38     filenames = d.GetFilenames()
39
40     # Get rid of any Warning while parsing the DICOM files
41     gdcm.Trace.WarningOff()
42
43     # instantiate Scanner:
44     sp = gdcm.Scanner.New();
45     s = sp.__ref__()
46     w = ProgressWatcher(s, 'Watcher')
47
48     s.AddTag( t1 );
49     s.AddTag( t2 );
50     b = s.Scan( filenames );
51     if(not b): sys.exit(1);
52
53     print "success" ;
54     #print s
55
56     pttv = gdcm.PythonTagToValue( s.GetMapping( filenames[1] ) )
57     pttv.Start()
58     # iterate until the end:
59     while( not pttv.IsAtEnd() ):
60         # get current value for tag and associated value:
61         # if tag was not found, then it was simply not added to the internal std::map
62         # Warning value can be None
63         tag = pttv.GetCurrentTag()
64         value = pttv.GetCurrentValue()
65         print tag,"->",value
66         # increment iterator
67         pttv.Next()

```

```

68
69  sys.exit(0)

```

## 29.141 SendFileSCU.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Perso/gdcm-gcc/bin
 * $ mono bin/SendFileSCU.exe server port input.dcm
 */
using System;
using gdcm;

public class SendFileSCU
{
    public static int Main(string[] args)
    {
        string server = args[0];
        ushort port = ushort.Parse(args[1]);
        string filename = args[2];

        bool b = CompositeNetworkFunctions.CEcho( server, port );
        if( !b ) return 1;

        FilenamesType files = new FilenamesType();
        files.Add( filename );
        b = CompositeNetworkFunctions.CStore( server, port, files );
        if( !b ) return 1;

        return 0;
    }
}

```

## 29.142 SimplePrint.cs

This is a C# example on how to use gdcm::SWIGDataSet

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
    Converter converter = new Converter();
    int a = converter.Convert<int>( some_int_blob );
    double b = converter.Convert<double>( some_double_blob );
 */

/*
 * Usage:

```

```

* $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
* $ mono bin/SimplePrint.exe gdcmData/012345.002.050.dcm
*/
using System;
using gdcm;

public class SimplePrint
{
    public static void RecurseDataSet(File f, DataSet ds, string indent)
    {
        CSharpDataSet cds = new CSharpDataSet(ds);
        while(!cds.IsAtEnd())
        {
            DataElement de = cds.GetCurrent();
            // Compute VR from the toplevel file, and the currently processed dataset:
            VR vr = DataSetHelper.ComputeVR(f, ds, de.GetTag() );

            if( vr.Compatible( new VR(VR.VRType.SQ) ) )
            {
                uint uvl = (uint)de.GetVL(); // Test cast is ok
                System.Console.WriteLine( indent + de.GetTag().toString() + ":" + uvl ); // why not ?
                //SequenceOfItems sq = de.GetSequenceOfItems();
                // GetValueAsSQ handle more cases than GetSequenceOfItems
                SmartPtrSQ sq = de.GetValueAsSQ();
                uint n = sq.GetNumberOfItems();
                for( uint i = 1; i <= n; i++) // item starts at 1, not 0
                {
                    Item item = sq.GetItem( i );
                    DataSet nested = item.GetNestedDataSet();
                    RecurseDataSet( f, nested, indent + " " );
                }
            }
            else
            {
                System.Console.WriteLine( indent + de.toString() );
            }
            cds.Next();
        }
    }

    public static int Main(string[] args)
    {
        string filename = args[0];
        Reader reader = new Reader();
        reader.SetFileName( filename );
        bool ret = reader.Read();
        if( !ret )
        {
            return 1;
        }
        File f = reader.GetFile();
        DataSet ds = f.GetDataSet();

        RecurseDataSet( f, ds, "" );

        return 0;
    }
}

```

## 29.143 SimplePrintPatientName.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
* Usage:
* $ export LD_LIBRARY_PATH=$HOME/Perso/gdcm/debug-gcc/bin

```

```

* $ mono bin/SimplePrintPatientName.exe gdcmlData/012345.002.050.dcm
*/
/*
This example was provided by Jonathan Morra /jonmorra gmail com/
on the gdcml mailing list (Fri, 28 May 2010)
*/
using System;
using gdcml;

namespace GDCMTest
{
    class SimplePrintPatientName
    {
        static int Main(string[] args)
        {
            if (args.Length != 1)
            {
                Console.WriteLine("This program prints the patient name of a dicom file with gdcml");
                Console.WriteLine("Usage: [input.dcm]");
                return 1;
            }

            gdcml.Reader reader = new gdcml.Reader();
            reader.SetFileName(args[0]);
            bool ret = reader.Read();
            //TagSetType tst = new TagSetType();
            //tst.Add( new Tag(0x7fe0,0x10) );
            //bool ret = reader.ReadUpToTag( new Tag(0x88,0x200), tst );
            if( !ret )
            {
                return 1;
            }

            gdcml.File file = reader.GetFile();

            gdcml.StringFilter filter = new gdcml.StringFilter();
            filter.SetFile(file);
            string value = filter.ToString(new gdcml.Tag(0x0010, 0x0010));

            Console.WriteLine("Patient Name: " + value);
            return 0;
        }
    }
}

```

## 29.144 SimpleScanner.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcml.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
* Simple example to show how to use Scanner API.
* It exposes the three different cases:
* - DICOM Attribute is present and has a value
* - DICOM Attribute is present and has no value
* - DICOM Attribute is not present at all
* It also shows the purpose of the function 'IsKey' to detect whether or
* not the file has been read by the gdcml::Scanner. Technically most of the time
* if a file is not a 'Key' this is because it is not a DICOM file. You need to use
* gdcml::System::FileExists to decide whether or not the file actually exist on the disk.
*
* It was tested on this particular image:
* ./SimpleScanner gdcmlData/012345.002.050.dcm
*/

#include "gdcmlScanner.h"
#include "gdcmlSimpleSubjectWatcher.h"

```

```

#include "gdcmFileNameEvent.h"

class MyFileWatcher : public gdcm::SimpleSubjectWatcher
{
public:
    MyFileWatcher(gdcm::Subject * s, const char *comment = ""):
        gdcm::SimpleSubjectWatcher(s,comment){}
    void ShowFileName(gdcm::Subject *, const gdcm::Event &evt)
    {
        const gdcm::FileNameEvent &pe = dynamic_cast<const
            gdcm::FileNameEvent&>(evt);
        const char *fn = pe.GetFileName();
        std::cout << "FileName: " << fn << " FileSize: " << gdcm::System::FileSize( fn )
            << std::endl;
    }
};

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        return 1;
    }
    const char *filename = argv[1];
    const char filename_invalid[] = "this is a file that may not exist on this disk.dcm";

    gdcm::SmartPointer<gdcm::Scanner> sp = new
        gdcm::Scanner;
    gdcm::Scanner &s = *sp;
    //gdcm::SimpleSubjectWatcher w(&s, "TestFileName" );
    MyFileWatcher w(&s, "TestFileName" );

    const gdcm::Tag tag_array[] = {
        gdcm::Tag(0x8,0x50),
        gdcm::Tag(0x8,0x51),
        gdcm::Tag(0x8,0x60),
    };
    s.AddTag( tag_array[0] );
    s.AddTag( tag_array[1] );
    s.AddTag( tag_array[2] );

    gdcm::Directory::FileNamesType filenames;
    filenames.push_back( filename );
    filenames.push_back( filename_invalid );

    if( !s.Scan( filenames ) )
    {
        return 1;
    }

    //s.Print( std::cout );

    if( s.IsKey( filename ) )
    {
        std::cout << "INFO:" << filename << " is a proper Key for the Scanner (this is a DICOM file)" <<
            std::endl;
    }

    if( !s.IsKey( filename_invalid ) )
    {
        std::cout << "INFO:" << filename_invalid << " is not a proper Key for the Scanner (this is either not a
            DICOM file or file does not exist)" << std::endl;
    }

    gdcm::Scanner::TagToValue const &ttv = s.GetMapping(filename);

    const gdcm::Tag *ptag = tag_array;
    for( ; ptag != tag_array + 3; ++ptag )
    {
        gdcm::Scanner::TagToValue::const_iterator it = ttv.find( *ptag );
        if( it != ttv.end() )
        {
            std::cout << *ptag << " was properly found in this file" << std::endl;
            // it contains a pair of value. the first one is the actual tag, so the following is always true:
            // *ptag == it->first
            // The second part is the actual value (stored as RAW strings). You will have to reinterpret this
            string
            // if VR for *ptag is not VR:VRASCII !
            const char *value = it->second;

```



```

        if( *value )
        {
            std::cout << "  It has the value: " << value << std::endl;
        }
        else
        {
            std::cout << "  It has no value (empty)" << std::endl;
        }
    }
    else
    {
        std::cout << "Sorry " << *ptag << " could not be found in this file" << std::endl;
    }
}

return 0;
}

```

## 29.145 SortImage.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
*/
#include "gdcmSorter.h"
#include "gdcmScanner.h"
#include "gdcmDataSet.h"
#include "gdcmAttribute.h"

bool mysort(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    //gdcm::Attribute<0x0020,0x0013> at1; // Instance Number
    gdcm::Attribute<0x0018,0x1060> at1; // Trigger Time
    gdcm::Attribute<0x0020,0x0032> at11; // Image Position (Patient)
    at1.Set( ds1 );
    at11.Set( ds1 );
    //gdcm::Attribute<0x0020,0x0013> at2;
    gdcm::Attribute<0x0018,0x1060> at2;
    gdcm::Attribute<0x0020,0x0032> at22;
    at2.Set( ds2 );
    at22.Set( ds2 );
    if( at11 == at22 )
    {
        return at1 < at2;
    }
    return at11 < at22;
}

bool mysort_part1(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    gdcm::Attribute<0x0018,0x1060> at1;
    at1.Set( ds1 );
    gdcm::Attribute<0x0018,0x1060> at2;
    at2.Set( ds2 );
    return at1 < at2;
}

bool mysort_part2(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    gdcm::Attribute<0x0020,0x0032> at1;
    at1.Set( ds1 );
    gdcm::Attribute<0x0020,0x0032> at2;
    at2.Set( ds2 );
    return at1 < at2;
}

```

```

// technically all files are in the same Frame of Reference, so this function
// should be a no-op
bool mysort_dummy(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    gdcm::Attribute<0x0020,0x0052> at1; // FrameOfReferenceUID
    at1.Set( ds1 );
    gdcm::Attribute<0x0020,0x0052> at2;
    at2.Set( ds2 );
    return at1 < at2;
}

int main(int argc, char *argv[])
{
    if (argc < 2 ) return 1;
    const char *dirname = argv[1];
    gdcm::Directory dir;
    unsigned int nfiles = dir.Load( dirname );

    dir.Print( std::cout );

    gdcm::Sorter sorter;
    sorter.SetSortFunction( mysort );
    sorter.Sort( dir.GetFilesNames() );

    std::cout << "Sorter:" << std::endl;
    sorter.Print( std::cout );

    gdcm::Sorter sorter2;
    sorter2.SetSortFunction( mysort_part1 );
    sorter2.StableSort( dir.GetFilesNames() );
    sorter2.SetSortFunction( mysort_part2 );
    sorter2.StableSort( sorter2.GetFilesNames() ); // IMPORTANT
    sorter2.SetSortFunction( mysort_dummy );
    sorter2.StableSort( sorter2.GetFilesNames() ); // IMPORTANT

    std::cout << "Sorter2:" << std::endl;
    sorter2.Print( std::cout );

    gdcm::Scanner s;
    s.AddTag( gdcm::Tag(0x20,0x32) ); // Image Position (Patient)
    //s.AddTag( gdcm::Tag(0x20,0x37) ); // Image Orientation (Patient)
    s.Scan( dir.GetFilesNames() );

    //s.Print( std::cout );

    // Count how many different IPP there are:
    const gdcm::Scanner::ValueType &values = s.GetValues();
    size_t nvalues = values.size();
    std::cout << "There are " << nvalues << " different type of values" << std::endl;

    //std::cout << "nfiles=" << nfiles << std::endl;
    if( nfiles % nvalues != 0 )
    {
        std::cerr << "Impossible: this is a not a proper series" << std::endl;
        return 1;
    }
    std::cout << "Series is composed of " << (nfiles/nvalues) << " different 3D volumes" << std::endl;

    return 0;
}

```

## 29.146 SortImage.py

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####

```

```

14
15 """
16 Usage:
17
18 python SortImage.py dirname
19 """
20
21 import gdc
22 import sys
23
24 def PrintProgress(object, event):
25     assert event == "ProgressEvent"
26     print "Progress:", object.GetProgress()
27
28 def MySort(ds1, ds2):
29     # compare ds1
30     return False
31
32 if __name__ == "__main__":
33
34     dirname = sys.argv[1]
35     d = gdc.Directory()
36     d.Load( dirname )
37
38     print d
39
40     sorter = gdc.Sorter()
41     sorter.SetSortFunction( MySort )
42     #sorter.AddObserver( "ProgressEvent", PrintProgress )
43     sorter.Sort( d.GetFilenames() )
44
45     print "Sorter:"
46     print sorter

```

## 29.147 SortImage2.cs

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdc.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/

/*
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdc/debug-gcc/bin
 * $ mono bin/SortImage.exe gdcData/012345.002.050.dcm out.dcm
 */
using System;
using gdc;

public class SortImage2
{
    bool mysort(DataSet ds1, DataSet ds2)
    {
        return false;
    }

    public static int Main(string[] args)
    {
        Sorter sorter = new Sorter();
        sorter.SetSortFunction( mysort );

        return 0;
    }
}

```

## 29.148 StandardizeFiles.cs

This is a C++ example on how to use ImageChangeTransferSyntax

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/

/*
 * Simple C# example to show how one would 'Standardize' a DICOM File-Set
 *
 * Usage:
 * $ export LD_LIBRARY_PATH=$HOME/Projects/gdcm/debug-gcc/bin
 * $ mono bin/StandardizeFiles.exe input_path output_path
 */
using System;
using gdcm;

public class StandardizeFiles
{
    public static bool ProcessOneFile( string filename, string outfilename )
    {
        PixmapReader reader = new PixmapReader();
        reader.SetFileName( filename );
        if( !reader.Read() )
        {
            System.Console.WriteLine( "Could not read: " + filename );
            return false;
        }

        ImageChangeTransferSyntax change = new ImageChangeTransferSyntax();
        change.SetForce( false ); // do we really want to recompress when input is already compressed in same
        // alg ?
        change.SetCompressIconImage( false ); // Keep it simple
        change.SetTransferSyntax( new TransferSyntax( TransferSyntax.TSType.JPEG2000Lossless ) );
        change.SetInput( reader.GetPixmap() );
        if( !change.Change() )
        {
            System.Console.WriteLine( "Could not change: " + filename );
            return false;
        }

        gdcm.FileMetaInformation fmi = reader.GetFile().GetHeader();
        // The following three lines make sure to regenerate any value:
        fmi.Remove( new gdcm.Tag(0x0002,0x0012) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0013) );
        fmi.Remove( new gdcm.Tag(0x0002,0x0016) );

        PixmapWriter writer = new PixmapWriter();
        writer.SetFileName( outfilename );
        writer.SetFile( reader.GetFile() );
        gdcm.Pixmap pixout = ((PixmapToPixmapFilter)change).GetOutput();

        writer.SetPixmap( pixout );
        if( !writer.Write() )
        {
            System.Console.WriteLine( "Could not write: " + outfilename );
            return false;
        }

        return true;
    }

    public static int Main(string[] args)
    {
        gdcm.FileMetaInformation.
            SetSourceApplicationEntityTitle( "My Standardize App" );

        // http://www.oid-info.com/get/1.3.6.1.4.17434
    }
}

```

```

string THERALYS_ORG_ROOT = "1.3.6.1.4.17434";
gdcm.UIDGenerator.SetRoot( THERALYS_ORG_ROOT );
System.Console.WriteLine( "Root dir is now: " + gdcm.UIDGenerator.
    GetRoot() );

string dir1 = args[0];
string dir2 = args[1];

// Check input is valid:
if( !gdcm.PosixEmulation.FileIsDirectory(dir1) )
{
    System.Console.WriteLine( "Input directory: " + dir1 + " does not exist. Sorry" );
    return 1;
}
if( !gdcm.PosixEmulation.FileIsDirectory(dir2) )
{
    System.Console.WriteLine( "Output directory: " + dir2 + " does not exist. Sorry" );
    return 1;
}

Directory d = new Directory();
uint nfiles = d.Load( dir1, true );
if(nfiles == 0) return 1;

// Process all filenames:
FilenameType filenames = d.GetFilesNames();
for( uint i = 0; i < nfiles; ++i )
{
    string filename = filenames[ (int)i ];
    string outfilename = filename.Replace( dir1, dir2 );
    System.Console.WriteLine( "Filename: " + filename );
    System.Console.WriteLine( "Out Filename: " + outfilename );
    if( !ProcessOneFile( filename, outfilename ) )
    {
        System.Console.WriteLine( "Could not process filename: " + filename );
        //return 1;
    }
}

return 0;
}

```

## 29.149 StreamImageReaderTest.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.
=====*/
// This work was realised during the GSOC 2011 by Manoj Alwani

#include "gdcmStreamImageReader.h"
#include "gdcmFileMetaInformation.h"
#include "gdcmSystem.h"
#include "gdcmFilename.h"
#include "gdcmByteSwap.h"
#include "gdcmTrace.h"
#include "gdcmTesting.h"
#include "gdcmImageHelper.h"
#include "gdcmImageReader.h"
#include "gdcmImage.h"
#include "gdcmMediaStorage.h"
#include "gdcmRAWCodec.h"
#include "gdcmJPEGLSCodec.h"
#include "gdcmUIDGenerator.h"
#include "gdcmStreamImageWriter.h"
#include "gdcmAttribute.h"

```

```

#include "gdcmFile.h"
#include "gdcmTag.h"

bool StreamImageRead(gdcm::StreamImageWriter & theStreamWriter,
    const char* filename, const char* outfilename, int resolution)
{
    gdcm::StreamImageReader reader;

    reader.SetFileName( filename );

    if (!reader.ReadImageInformation())
    {
        std::cerr << "unable to read image information" << std::endl;
        return 1; //unable to read tags as expected.
    }
    //let's be tricky; each image will be read in portions, first the top half, then the bottom
    //that way, we can test how the stream handles fragmentation of the data
    //we could also loop this to get various different size combinations, but I'm not sure
    //that's useful, yet.
    std::vector<unsigned int> extent =
        gdcm::ImageHelper::GetDimensionsValue(reader.
            GetFile());
    // std::cout << extent[0];
    //at this point, these values aren't used, but may be in the future
    //unsigned short xmin = 0;
    //unsigned short xmax = extent[0];
    //unsigned short ymin = 0;
    //unsigned short ymax = extent[1];
    //unsigned short zmin = 0;
    //unsigned short zmax = extent[2];

    std::cout<< "\n Row: "<<extent[0] <<"\n Col :"<< extent[1]<< "\n Resolution :"<< extent[2] << std::endl;

    int a =1;
    for (int i=1; i<=(extent[2]-resolution);++i)
        a = a*2;

    reader.DefinePixelExtent(0, extent[0]/a, 0, extent[1]/a, resolution-1, resolution);

    unsigned long len = reader.DefineProperBufferLength();
    char* finalBuffer = new char[len];
    memset(finalBuffer, 0, sizeof(char)*len);

    if (reader.CanReadImage())
    {
        bool result = reader.Read(finalBuffer, len);
        if( !result )
        {
            std::cout << "res2 failure:" << filename << std::endl;
            delete [] finalBuffer;
            return 1;
        }
        else
        {
            std::cout<< "Able to read";
        }
    }
    else
    {
        std::cerr<< "Not able to put in buffer"<< std::endl;
    }
}

/*
    //now, read in smaller buffer extents
    reader.DefinePixelExtent(xmin, xmax, ymin, ymax);
    len = reader.DefineProperBufferLength();

    char* buffer = new char[len];
    bool res2 = reader.Read(buffer, len);
    if( !res2 ){
        std::cerr << "res2 failure:" << filename << std::endl;
        return 1;
    }
    //copy the result into finalBuffer
    memcpy(finalBuffer, buffer, len);

    //now read the next half of the image
    ymin = ymax;
    ymax = extent[1];

    reader.DefinePixelExtent(xmin, xmax, ymin, ymax);

```

```

//std::cerr << "Success to read image from file: " << filename << std::endl;
unsigned long len2 = reader.DefineProperBufferLength();

char* buffer2 = new char[len2];
bool res3 = reader.Read(buffer2, len2);
if( !res3 ){
    std::cerr << "res3 failure:" << filename << std::endl;
    return 1;
}
//copy the result into finalBuffer
memcpy(&(finalBuffer[len]), buffer2, len2);

delete [] buffer;
delete [] buffer2;
*/

gdcm::Writer w;
gdcm::File &file = w.GetFile();
gdcm::DataSet &ds = file.GetDataSet();

file.GetHeader().SetDataSetTransferSyntax(
    gdcm::TransferSyntax::ExplicitVRLittleEndian );

gdcm::UIDGenerator uid;
gdcm::DataElement de( gdcm::Tag(0x8,0x18) ); // SOP Instance UID
de.SetVR( gdcm::VR::UI );
const char *u = uid.Generate();
de.SetByteValue( u, strlen(u) );
ds.Insert( de );

gdcm::DataElement del( gdcm::Tag(0x8,0x16) );
del.SetVR( gdcm::VR::UI );
gdcm::MediaStorage ms(
    gdcm::MediaStorage::VLWholeSlideMicroscopyImageStorage
);
del.SetByteValue( ms.GetString(), strlen(ms.GetString()) );
ds.Insert( del );

const char mystr[] = "MONOCHROME2 ";
gdcm::DataElement de2( gdcm::Tag(0x28,0x04) );
//de.SetTag(gdcm::Tag(0x28,0x04));
de2.SetVR( gdcm::VR::CS );
de2.SetByteValue(mystr, strlen(mystr));
ds.Insert( de2 );

gdcm::Attribute<0x0028,0x0008> Number_Of_Frames = {1};
ds.Insert( Number_Of_Frames.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0010> row = {extent[0]/a}; //
ds.Insert( row.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0011> col = {extent[1]/a}; //
ds.Insert( col.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0100> at = {8};
ds.Insert( at.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0002> at1 = {1}; //
ds.Insert( at1.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0101> at2 = {8};
ds.Insert( at2.GetAsDataElement() );

gdcm::Attribute<0x0028,0x0102> at3 = {7};
ds.Insert( at3.GetAsDataElement() );
/*
ds1.Remove( gdcm::Tag(0x0028,0x0008) );

gdcm::Attribute<0x0028,0x0008> Number_Of_Frames = {1};
ds1.Insert( Number_Of_Frames.GetAsDataElement() );
*/
theStreamWriter.SetFile(file);

if (!theStreamWriter.WriteImageInformation())
{
    std::cerr << "unable to write image information" << std::endl;
    return 1; //the CanWrite function should prevent getting here, else,
    //that's a test failure
}
std::vector<unsigned int> extent1 = gdcm::ImageHelper::GetDimensionsValue
(file);

```

```

unsigned short xmax = extent1[0];
unsigned short ymax = extent1[1];
unsigned short theChunkSize = 1;
unsigned short ychunk = extent1[1]/theChunkSize; //go in chunk sizes of theChunkSize
unsigned short zmax = 1;

std::cout<< "\n Row: "<<extent1[0] <<"\n Col : "<< extent1[1]<< "\n Resolution : "<< extent1[2] <<
std::endl;

if (xmax == 0 || ymax == 0)
{
std::cerr << "Image has no size, unable to write zero-sized image." << std::endl;
return 0;
}

int z, y, nexty;
unsigned long prevLen = 0; //when going through the char buffer, make sure to grab
//the bytes sequentially. So, store how far you got in the buffer with each iteration.

for (z = 0; z < zmax; ++z){
for (y = 0; y < ymax; y += ychunk){
nexty = y + ychunk;
if (nexty > ymax) nexty = ymax;
theStreamWriter.DefinePixelExtent(0, xmax, y, nexty, z, z+1);
unsigned long len = theStreamWriter.DefineProperBufferLength();
std::cout << "\n" <<len;
char* finalBuffer1 = new char[len];
memcpy(finalBuffer1, &(finalBuffer[prevLen]), len);
std::cout << "\nable to write";

if (!theStreamWriter.Write(finalBuffer1, len)){
std::cerr << "writing failure:" << "output.dcm" << " at y = " << y << " and z= " << z <<
std::endl;
delete [] finalBuffer1;
delete [] finalBuffer;
return 1;
}
delete [] finalBuffer1;
prevLen += len;
}
}
delete [] finalBuffer;
std::cout << "all is set";

return true;
}

int main(int argc, char *argv[])
{
if( argc < 3 )
{
std::cerr << argv[0] << " input.dcm output.dcm Resolution" << std::endl;
return 1;
}

const char *filename = argv[1];
const char *outfilename = argv[2];
char *res = argv[3];

int resolution = atoi(res);

gdcm::StreamImageWriter theStreamWriter;

std::ofstream of;
of.open( outfile, std::ios::out | std::ios::binary );
theStreamWriter.SetStream(of);

// else
// First of get rid of warning/debug message
gdcm::Trace::DebugOn();
gdcm::Trace::WarningOn();

if(!StreamImageRead( theStreamWriter, filename, outfile, resolution))
return 1;

uint16_t firstTag1 = 0xfffe;
uint16_t secondTag1 = 0xe0dd;
uint32_t thirdTag1 = 0x00000000;
//uint16_t fourthTag1 = 0xffff;

```



```

const int theBufferSize1 = 2*sizeof(uint16_t)+sizeof(uint32_t);
char* tmpBuffer2 = new char[theBufferSize1];
memcpy(&(tmpBuffer2[0]), &firstTag1, sizeof(uint16_t));
memcpy(&(tmpBuffer2[sizeof(uint16_t)]), &secondTag1, sizeof(uint16_t));
memcpy(&(tmpBuffer2[2*sizeof(uint16_t)]), &thirdTag1, sizeof(uint32_t));
//memcpy(&(tmpBuffer2[3*sizeof(uint16_t)]), &fourthTag1, sizeof(uint16_t));
assert( of && !of.eof() && of.good() );
of.write(tmpBuffer2, theBufferSize1);
of.flush();
assert( of );

return 0;
}

```

## 29.150 TestByteSwap.cxx

This is a C++ example on how to use `gdcm::ByteSwap`

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmTypes.h"
#include "gdcmSwapCode.h"
#include "gdcmByteSwap.h"

#include <string.h> // memcpy

int myfunc()
{
    char vl_str[4];
    const char raw[] = "\000\000\000\004";
    memcpy(vl_str, raw, 4);
    uint32_t vl;
    gdcm::ByteSwap<uint32_t>::SwapRangeFromSwapCodeIntoSystem(
        ((uint32_t*)(&vl_str)), gdcm::SwapCode::BigEndian, 1);
    memcpy(&vl, vl_str, 4);
    if( vl != 0x00000004 )
    {
        std::cerr << std::hex << "vl: " << vl << std::endl;
        return 1;
    }

    gdcm::ByteSwap<uint32_t>::SwapFromSwapCodeIntoSystem(
        vl, gdcm::SwapCode::LittleEndian);
    if( vl != 0x00000004 )
    {
        std::cerr << std::hex << "vl: " << vl << std::endl;
        return 1;
    }

    gdcm::ByteSwap<uint32_t>::SwapFromSwapCodeIntoSystem(
        vl, gdcm::SwapCode::BigEndian);
    std::cout << std::hex << "vl: " << vl << std::endl;
    if( vl != 0x4000000 )
    {
        return 1;
    }

    return 0;
}

int TestByteSwap(int , char *[])
{
    gdcm::SwapCode sc = gdcm::SwapCode::Unknown;
    if ( gdcm::ByteSwap<uint16_t>::SystemIsBigEndian() )
    {

```

```

    sc = gdcmm::SwapCode::BigEndian;
}
else if ( gdcmm::ByteSwap<uint16_t>::SystemIsLittleEndian()
)
{
    sc = gdcmm::SwapCode::LittleEndian;
}
if( sc == gdcmm::SwapCode::Unknown )
{
    return 1;
}

std::cout << "sc: " << sc << std::endl;

uint16_t t = 0x1234;
gdcmm::ByteSwap<uint16_t>::SwapFromSwapCodeIntoSystem(
    t, sc);
if( sc == gdcmm::SwapCode::BigEndian )
{
    if( t != 0x3412 )
    {
        std::cerr << std::hex << "t: " << t << std::endl;
        return 1;
    }
    // ok test pass rest value to old one
    t = 0x1234;
}
else if ( sc == gdcmm::SwapCode::LittleEndian )
{
    if( t != 0x1234 )
    {
        std::cerr << std::hex << "t: " << t << std::endl;
        return 1;
    }
}

union { char n[2]; uint16_t tn; } ul6;
memcpy(ul6.n, &t, 2 );
gdcmm::ByteSwap<uint16_t>::SwapRangeFromSwapCodeIntoSystem
    (&ul6.tn, sc, 1);
uint16_t tn = ul6.tn;
if( sc == gdcmm::SwapCode::BigEndian )
{
    if( tn != 0x3412 )
    {
        std::cerr << std::hex << "tn: " << tn << std::endl;
        return 1;
    }
    // ok test pass rest value to old one
    t = 0x1234;
}
else if ( sc == gdcmm::SwapCode::LittleEndian )
{
    if( tn != 0x1234 )
    {
        std::cerr << std::hex << "tn: " << tn << std::endl;
        return 1;
    }
}
gdcmm::ByteSwap<uint16_t>::SwapRangeFromSwapCodeIntoSystem
    (&ul6.tn, gdcmm::SwapCode::BigEndian, 1);
tn = ul6.tn;
if( sc == gdcmm::SwapCode::LittleEndian )
{
    if( tn != 0x3412 )
    {
        std::cerr << std::hex << "tn: " << tn << std::endl;
        return 1;
    }
}
else if ( sc == gdcmm::SwapCode::BigEndian )
{
    if( tn != 0x1234 )
    {
        std::cerr << std::hex << "tn: " << tn << std::endl;
        return 1;
    }
}

if( myfunc() )
{

```

```

        return 1;
    }

    uint16_t array[] = { 0x1234 };
    gdcm::ByteSwap<uint16_t>::SwapRangeFromSwapCodeIntoSystem
        (array,
         gdcm::SwapCode::BigEndian,2);
    if ( array[0] != 0x3412 )
    {
        return 1;
    }

    return 0;
}

```

## 29.151 TestReader.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmReader.h"
#include "gdcmFileMetaInformation.h"
#include "gdcmFile.h"
#include "gdcmTesting.h"
#include "gdcmMediaStorage.h"

int TestRead(const char* filename, bool verbose = false)
{
    if( verbose )
        std::cout << "TestRead: " << filename << std::endl;

    gdcm::Reader reader;
    reader.SetFileName( filename );
    if ( !reader.Read() )
    {
        std::cerr << "TestReadError: Failed to read: " << filename << std::endl;
        return 1;
    }

    //commenting out the fmi and ds to avoid warnings
    //const gdcm::FileMetaInformation &h = reader.GetFile().GetHeader();
    //std::cout << h << std::endl;

    //const gdcm::DataSet &ds = reader.GetFile().GetDataSet();
    //std::cout << ds << std::endl;

    const char *ref = gdcm::Testing::GetMediaStorageFromFile(filename);
    gdcm::MediaStorage ms;
    ms.SetFromFile( reader.GetFile() );
    if( !ref )
    {
        std::cerr << "TestReadError: Missing MediaStorage: " << filename << std::endl;
        std::cerr << "It should be: " << ms << std::endl;
        return 1;
    }

    if( ms.IsUndefined() && ref && *ref != 0 )
    {
        std::cerr << "TestReadError: MediaStorage: " << filename << std::endl;
        std::cerr << "It should be instead: " << ref << std::endl;
        return 1;
    }

    // Make sure it is the right one:

    if( ref && *ref != 0 && ms != gdcm::MediaStorage::GetMSType(ref) )
    {

```

```

        std::cerr << "Error: Found MediaStorage: " << ms << " for " << filename << std::endl;
        std::cerr << "It should be instead: " << ref << std::endl;
        return 1;
    }

    return 0;
}

int TestReader(int argc, char *argv[])
{
    if( argc == 2 )
    {
        const char *filename = argv[1];
        return TestRead(filename, true);
    }

    // else
    gdc::Trace::DebugOff();
    gdc::Trace::WarningOff();
    int r = 0, i = 0;
    const char *filename;
    const char * const *filenames = gdc::Testing::GetFileNames();
    while( (filename = filenames[i]) )
    {
        r += TestRead( filename );
        ++i;
    }

    return r;
}

```

## 29.152 TestReader.py

This is a C++ example on how to use [gdc::Reader](#)

```

1 #####
2 #
3 #   Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 #   Copyright (c) 2006-2011 Mathieu Malaterre
6 #   All rights reserved.
7 #   See Copyright.txt or http://gdc.sourceforge.net/Copyright.html for details.
8 #
9 #   This software is distributed WITHOUT ANY WARRANTY; without even
10 #   the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 #   PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 import gdc
16 import os,sys
17
18 def TestRead(filename, verbose = False):
19     r = gdc.Reader()
20     r.SetFileName( filename )
21     sucess = r.Read()
22     #if verbose: print r.GetFile()
23     if verbose: print(r.GetFile().GetDataSet())
24     return sucess
25
26 if __name__ == "__main__":
27     sucess = 0
28     try:
29         filename = os.sys.argv[1]
30         sucess += TestRead( filename, True )
31     except:
32         # loop over all files:
33         gdc.Trace.DebugOff()
34         gdc.Trace.WarningOff()
35         t = gdc.Testing()
36         nfiles = t.GetNumberOfFileNames()
37         for i in range(0,nfiles):
38             filename = t.GetFileName(i)
39             sucess += TestRead( filename )
40
41

```

```

42  # Test succeed ?
43  sys.exit(sucess == 0)

```

## 29.153 threadgdcmm.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmmReader.h"
#include "gdcmmImageReader.h"
#include "gdcmmDirectory.h"
#include "gdcmmSystem.h"

#include "vtkImageData.h"
#include "vtkStructuredPointsWriter.h"

#include <pthread.h>

struct threadparams
{
    const char **filenames;
    size_t nfiles;
    char *scalarpointer;
// TODO I should also pass in the dim of the reference image just in case
};

void *ReadFilesThread(void *voidparams)
{
    const threadparams *params = static_cast<const threadparams *> (voidparams);

    const size_t nfiles = params->nfiles;
    for(unsigned int file = 0; file < nfiles; ++file)
    {
        /*
        // TODO: update progress
        pthread_mutex_lock(&params->lock);
        //section critique
        ReadingProgress+=params->stepProgress;
        pthread_mutex_unlock(&params->lock);
        */
        const char *filename = params->filenames[file];
        //std::cerr << filename << std::endl;

        gdcmm::ImageReader reader;
        reader.SetFileName( filename );
        try
        {
            {
                if( !reader.Read() )
                {
                    std::cerr << "Failed to read: " << filename << std::endl;
                    break;
                }
            }
        } catch( ... )
        {
            std::cerr << "Failed to read: " << filename << std::endl;
            break;
        }

        const gdcmm::Image &image = reader.GetImage();
        unsigned long len = image.GetBufferLength();
        char * pointer = params->scalarpointer;
        #if 0
        char *tempimage = new char[len];
        image.GetBuffer(tempimage);

        memcpy(pointer + file*len, tempimage, len);

```

```

        delete[] tempimage;
    #else
        char *tempimage = pointer + file * len;
        image.GetBuffer(tempimage);
    #endif
    }

    return voidparams;
}

void ShowFileNames(const threadparams &params)
{
    std::cout << "start" << std::endl;
    for(unsigned int i = 0; i < params.nfiles; ++i)
    {
        const char *filename = params filenames[i];
        std::cout << filename << std::endl;
    }
    std::cout << "end" << std::endl;
}

void ReadFiles(size_t nfiles, const char *filenames[])
{
    // \precondition: nfiles > 0
    assert( nfiles > 0 );
    const char *reference= filenames[0]; // take the first image as reference

    gdcm::ImageReader reader;
    reader.SetFileName( reference );
    if( !reader.Read() )
    {
        // That would be very bad...
        assert(0);
    }

    const gdcm::Image &image = reader.GetImage();
    gdcm::PixelFormat pixeltype = image.GetPixelFormat();
    unsigned long len = image.GetBufferLength();
    const unsigned int *dims = image.GetDimensions();
    unsigned short pixelsize = pixeltype.GetPixelSize();
    (void)pixelsize;
    assert( image.GetNumberOfDimensions() == 2 );

    vtkImageData *output = vtkImageData::New();
    output->SetDimensions(dims[0], dims[1], (int)nfiles);

    #if (VTK_MAJOR_VERSION >= 6)
    int numscal = pixeltype.GetSamplesPerPixel();
    switch( pixeltype )
    {
        case gdcm::PixelFormat::INT8:
            output->AllocateScalars( VTK_SIGNED_CHAR, numscal );
            break;
        case gdcm::PixelFormat::UINT8:
            output->AllocateScalars( VTK_UNSIGNED_CHAR, numscal );
            break;
        case gdcm::PixelFormat::INT16:
            output->AllocateScalars( VTK_SHORT, numscal );
            break;
        case gdcm::PixelFormat::UINT16:
            output->AllocateScalars( VTK_UNSIGNED_SHORT, numscal );
            break;
        case gdcm::PixelFormat::INT32:
            output->AllocateScalars( VTK_INT, numscal );
            break;
        case gdcm::PixelFormat::UINT32:
            output->AllocateScalars( VTK_UNSIGNED_INT, numscal );
            break;
        default:
            assert(0);
    }
    #else
    switch( pixeltype )
    {
        case gdcm::PixelFormat::INT8:
            #if (VTK_MAJOR_VERSION >= 5) || ( VTK_MAJOR_VERSION == 4 && VTK_MINOR_VERSION > 5 )
            output->SetScalarType ( VTK_SIGNED_CHAR );
            #else
            output->SetScalarType ( VTK_CHAR );
            #endif
            break;
    }
    #endif
}

```

```

case gdcm::PixelFormat::UINT8:
    output->SetScalarType ( VTK_UNSIGNED_CHAR );
    break;
case gdcm::PixelFormat::INT16:
    output->SetScalarType ( VTK_SHORT );
    break;
case gdcm::PixelFormat::UINT16:
    output->SetScalarType ( VTK_UNSIGNED_SHORT );
    break;
case gdcm::PixelFormat::INT32:
    output->SetScalarType ( VTK_INT );
    break;
case gdcm::PixelFormat::UINT32:
    output->SetScalarType ( VTK_UNSIGNED_INT );
    break;
default:
    assert(0);
}
output->SetNumberOfScalarComponents ( pixeltype.GetSamplesPerPixel() );
output->AllocateScalars();
#endif
char * scalarpointer = static_cast<char*>(output->GetScalarPointer());

const unsigned int nthreads = 4;
threadparams params[nthreads];

//pthread_mutex_t lock;
//pthread_mutex_init(&lock, NULL);

pthread_t *pthread = new pthread_t[nthreads];

// There is nfiles, and nThreads
assert( nfiles > nthreads );
const size_t partition = nfiles / nthreads;
for (unsigned int thread=0; thread < nthreads; ++thread)
{
    params[thread].filenames = filenames + thread * partition;
    params[thread].nfiles = partition;
    if( thread == nthreads - 1 )
    {
        // There is slightly more files to process in this thread:
        params[thread].nfiles += nfiles % nthreads;
    }
    assert( thread * partition < nfiles );
    params[thread].scalarpointer = scalarpointer + thread * partition * len;
    //assert( params[thread].scalarpointer < scalarpointer + 2 * dims[0] * dims[1] * dims[2] );
    // start thread:
    int res = pthread_create( &pthread[thread], NULL, ReadFilesThread, &params[thread]);
    if( res )
    {
        std::cerr << "Unable to start a new thread, pthread returned: " << res << std::endl;
        assert(0);
    }
    //ShowFileNames(params[thread]);
}
// DEBUG
size_t total = 0;
for (unsigned int thread=0; thread < nthreads; ++thread)
{
    total += params[thread].nfiles;
}
assert( total == nfiles );
// END DEBUG

for (unsigned int thread=0; thread<nthreads;thread++)
{
    pthread_join( pthread[thread], NULL);
}
delete[] pthread;

//pthread_mutex_destroy(&lock);

// For some reason writing down the file is painfully slow...
vtkStructuredPointsWriter *writer = vtkStructuredPointsWriter::New();
#if (VTK_MAJOR_VERSION >= 6)
    writer->SetInputData( output );
#else
    writer->SetInput( output );
#endif
writer->SetFileName( "/tmp/threadgdcm.vtk" );
writer->SetFileTypeToBinary();

```

```

//writer->Write();
writer->Delete();

//output->Print( std::cout );
output->Delete();
}

int main(int argc, char *argv[])
{
    if( argc < 2 )
    {
        std::cerr << argv[0] << " [directory|list of filenames]\n";
        return 1;
    }

    // Check if user pass in a single directory
    if( argc == 2 && gdcm::System::FileIsDirectory( argv[1] ) )
    {
        gdcm::Directory d;
        d.Load( argv[1] );
        gdcm::Directory::FileNamesType l = d.
            GetFileNames();
        const size_t nfiles = l.size();
        const char **filenames = new const char* [ nfiles ];
        for(unsigned int i = 0; i < nfiles; ++i)
        {
            filenames[i] = l[i].c_str();
        }
        ReadFiles(nfiles, filenames);
        delete[] filenames;
    }
    else
    {
        // Simply copy all filenames into the vector:
        const char **filenames = const_cast<const char**>(argv+1);
        const size_t nfiles = argc - 1;
        ReadFiles(nfiles, filenames);
    }

    return 0;
}

```

## 29.154 TraverseModules.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
*/

#include "gdcmDefs.h"
#include "gdcmGlobal.h"
#include "gdcmIODs.h"
#include "gdcmIOD.h"
#include "gdcmMacros.h"
#include "gdcmIODEntry.h"
#include "gdcmModules.h"
#include "gdcmModule.h"
#include "gdcmAnonymizer.h"
#include "gdcmDicts.h"

int main(int , char *[])
{
    using namespace gdcm;
    static Global &g = Global::GetInstance();

```



```

if( !g.LoadResourcesFiles() )
{
    return 1;
}

static const Defs &defs = g.GetDefs();
static const Modules &modules = defs.GetModules();
static const IODs &iods = defs.GetIODs();
static const Macros &macros = defs.GetMacros();
static const Dicts &dicts = g.GetDicts();

std::vector<Tag> tags =
    gdcmm::Anonymizer::GetBasicApplicationLevelConfidentialityProfileAttributes
    ();
for( std::vector<Tag>::const_iterator tit = tags.begin(); tit != tags.end(); ++tit )
{
    const Tag &tag = *tit;
    const DictEntry &dictentry = dicts.GetDictEntry(tag);
    std::cout << "Processing Attribute: " << tag << " " << dictentry << std::endl;

    IODs::IODMapTypeConstIterator it = iods.Begin();
    for( ; it != iods.End(); ++it )
    {
        const IODs::IODName &name = it->first;
        const IOD &iod = it->second;

        const size_t niods = iod.GetNumberOfIODs();
        // Iterate over each iod entry in order:
        for(unsigned int idx = 0; idx < niods; ++idx)
        {
            const IODEntry &iodentry = iod.GetIODEntry(idx);
            const char *ref = iodentry.GetRef();
            //Usage::UsageType ut = iodentry.GetUsageType();

            const Module &module = modules.GetModule( ref );
            if( module.FindModuleEntryInMacros(macros, tag) )
            {
                const ModuleEntry &module_entry = module.
                    GetModuleEntryInMacros(macros,tag);
                Type type = module_entry.GetType();
                std::cout << "IOD Name: " << name << std::endl;
                std::cout << "Type: " << type << std::endl;
            }
        }
    }
}

return 0;
}

```

## 29.155 uid\_unique.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcmm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
#include "gdcmmUIDGenerator.h"

#include <iostream>
#include <string>
#include <set>

int main()
{
    gdcmm::UIDGenerator uid;
    //const char myroot[] = "9876543210.9876543210.9876543210.9876543210.9876543210"; // fails in ~40000
    tries

```

```

const char myroot[] = "9876543210.9876543210.9876543210";
uid.SetRoot( myroot );
std::set<std::string> uids;
uint64_t wrap = 0;
uint64_t c = 0;
while(1)
{
    const char *unique = uid.Generate();
    //std::cout << unique << std::endl;
    if( c % 10000 == 0 )
    {
        std::cout << "wrap=" << wrap << ",c=" << c << std::endl;
    }
    ++c;
    if( c == 0 )
    {
        wrap++;
    }
    if ( uids.count(unique) == 1 )
    {
        std::cerr << "Failed with: " << unique << std::endl;
        return 1;
    }
    uids.insert( unique );
}
}

```

## 29.156 VolumeSorter.cxx

```

/*=====
Program: GDCM (Grassroots DICOM). A DICOM library

Copyright (c) 2006-2011 Mathieu Malaterre
All rights reserved.
See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.

This software is distributed WITHOUT ANY WARRANTY; without even
the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
PURPOSE. See the above copyright notice for more information.

=====*/
/*
*/
#include "gdcmSorter.h"
#include "gdcmIPPSorter.h"
#include "gdcmScanner.h"
#include "gdcmDataSet.h"
#include "gdcmAttribute.h"
#include "gdcmTesting.h"

bool mysort1(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    gdcm::Attribute<0x0020,0x000d> at1;
    at1.Set( ds1 );
    gdcm::Attribute<0x0020,0x000d> at2;
    at2.Set( ds2 );
    return at1 < at2;
}

bool mysort2(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    gdcm::Attribute<0x0020,0x000e> at1;
    at1.Set( ds1 );
    gdcm::Attribute<0x0020,0x000e> at2;
    at2.Set( ds2 );
    return at1 < at2;
}

bool mysort3(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    // This is a floating point number is the comparison ok ?
    gdcm::Attribute<0x0020,0x0037> at1;
    at1.Set( ds1 );
    gdcm::Attribute<0x0020,0x0037> at2;
    at2.Set( ds2 );
}

```

```

    return at1 < at2;
}

bool mysort4(gdcm::DataSet const & ds1, gdcm::DataSet const & ds2 )
{
    // Do the IPP sorting here
    gdcm::Attribute<0x0020,0x0032> iopl1;
    gdcm::Attribute<0x0020,0x0037> iopl;
    iopl1.Set( ds1 );
    iopl.Set( ds1 );
    gdcm::Attribute<0x0020,0x0032> iopl2;
    gdcm::Attribute<0x0020,0x0037> iopl2;
    iopl2.Set( ds2 );
    iopl2.Set( ds2 );
    if( iopl1 != iopl2 )
    {
        return false;
    }

    // else
    double normal[3];
    normal[0] = iopl1[1]*iopl1[5] - iopl1[2]*iopl1[4];
    normal[1] = iopl1[2]*iopl1[3] - iopl1[0]*iopl1[5];
    normal[2] = iopl1[0]*iopl1[4] - iopl1[1]*iopl1[3];
    double dist1 = 0;
    for( int i = 0; i < 3; ++i) dist1 += normal[i]*iopl1[i];
    double dist2 = 0;
    for( int i = 0; i < 3; ++i) dist2 += normal[i]*iopl2[i];

    std::cout << dist1 << ", " << dist2 << std::endl;
    return dist1 < dist2;
}

int main(int argc, char *argv[])
{
    const char *extradataroot = gdcm::Testing::GetDataExtraRoot();
    std::string dirl;
    if( argc < 2 )
    {
        if( !extradataroot )
        {
            return 1;
        }
        dirl = extradataroot;
        dirl += "/gdcmSampleData/ForSeriesTesting/VariousIncidences/ST1";
    }
    else
    {
        dirl = argv[1];
    }

    gdcm::Directory d;
    d.Load( dirl.c_str(), true ); // recursive !
    const gdcm::Directory::FileNamesType &ll = d.
        GetFileNames();
    const size_t nfiles = ll.size();
    std::cout << nfiles << std::endl;

    //if( nfiles != 280 )
    // {
    //     return 1;
    // }

    //d.Print( std::cout );

    gdcm::Scanner s0;
    const gdcm::Tag t1(0x0020,0x000d); // Study Instance UID
    const gdcm::Tag t2(0x0020,0x000e); // Series Instance UID
    //const gdcm::Tag t3(0x0010,0x0010); // Patient's Name
    s0.AddTag( t1 );
    s0.AddTag( t2 );
    //s0.AddTag( t3 );
    //s0.AddTag( t4 );
    //s0.AddTag( t5 );
    //s0.AddTag( t6 );
    bool b = s0.Scan( d.GetFileNames() );
    if( !b )
    {
        std::cerr << "Scanner failed" << std::endl;
    }
}

```

```

    return 1;
}

//s0.Print( std::cout );

// Only get the DICOM files:
gdcm::Directory::FileNamesType l2 = s0.GetKeys();
const size_t nfiles2 = l2.size();
std::cout << nfiles2 << std::endl;

if ( nfiles2 > nfiles )
{
    return 1;
}

gdcm::Sorter sorter;
sorter.SetSortFunction( mysort1 );
sorter.StableSort( l2 );

sorter.SetSortFunction( mysort2 );
sorter.StableSort( sorter.GetFileNames() );

sorter.SetSortFunction( mysort3 );
sorter.StableSort( sorter.GetFileNames() );

sorter.SetSortFunction( mysort4 );
sorter.StableSort( sorter.GetFileNames() );

//sorter.Print( std::cout );

// Let's try to check our result:
// assume that IPP is precise enough so that we can test floating point equality:
size_t nvalues = 0;
{
    gdcm::Scanner s;
    s.AddTag( gdcm::Tag(0x20,0x32) ); // Image Position (Patient)
    //s.AddTag( gdcm::Tag(0x20,0x37) ); // Image Orientation (Patient)
    s.Scan( d.GetFileNames() );

    //s.Print( std::cout );

    const gdcm::Scanner::ValuesType &values = s.GetValues();
    nvalues = values.size();
    std::cout << "There are " << nvalues << " different type of values" << std::endl;
    assert( nfiles2 % nvalues == 0 );
    std::cout << "Series is composed of " << (nfiles/nvalues) << " different 3D volumes" << std::endl;
}

gdcm::Directory::FileNamesType sorted_files = sorter.
    GetFileNames();

// Which means we can take nvalues files at a time and execute gdcm::IPPSorter on it:
gdcm::IPPSorter ippsorter;
gdcm::Directory::FileNamesType sub( sorted_files.begin(), sorted_files.
    begin() + nvalues);
std::cout << sub.size() << std::endl;
std::cout << sub[0] << std::endl;
std::cout << sub[nvalues-1] << std::endl;
ippsorter.SetComputeZSpacing( false );
if( !ippsorter.Sort( sub ) )
{
    std::cerr << "Could not sort" << std::endl;
    return 1;
}

std::cout << "IPPSorter:" << std::endl;
ippsorter.Print( std::cout );

return 0;
}

```

## 29.157 WriteBuffer.py

```

1 #####
2 #

```

```

3 # Program: GDCM (Grassroots DICOM). A DICOM library
4 #
5 # Copyright (c) 2006-2011 Mathieu Malaterre
6 # All rights reserved.
7 # See Copyright.txt or http://gdcm.sourceforge.net/Copyright.html for details.
8 #
9 # This software is distributed WITHOUT ANY WARRANTY; without even
10 # the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR
11 # PURPOSE. See the above copyright notice for more information.
12 #
13 #####
14
15 """
16 Usage:
17
18 http://chuckhahm.com/Ischem/Zurich/XX_0134
19
20 (2005,1132) SQ (Sequence with undefined length #=8) # u/1, 1 Unknown Tag & Data
21 (fffe,e000) na (Item with undefined length #=9) # u/1, 1 Item
22 (2005,0011) LO [Philips MR Imaging DD 002] # 26, 1 PrivateCreator
23 (2005,1137) PN [PDF_CONTROL_GEN_PARS] # 20, 1 Unknown Tag & Data
24 (2005,1138) PN (no value available) # 0, 0 Unknown Tag & Data
25 (2005,1139) PN [IEEE_PDF] # 8, 1 Unknown Tag & Data
26 (2005,1140) PN (no value available) # 0, 0 Unknown Tag & Data
27 (2005,1141) PN (no value available) # 0, 0 Unknown Tag & Data
28 (2005,1143) SL 3103 # 4, 1 Unknown Tag & Data
29 (2005,1144) OW 0566\0000\013b\0000\0a4a\0000\000e\0000\0a7a\0000\0195\0000\0008... # 3104, 1 Unknown
Tag & Data
30 (2005,1147) CS [Y] # 2, 1 Unknown Tag & Data
31 (fffe,e00d) na (ItemDelimitationItem) # 0, 0 ItemDelimitationItem
32 (fffe,e000) na (Item with undefined length #=9) # u/1, 1 Item
33 (2005,0011) LO [Philips MR Imaging DD 002] # 26, 1 PrivateCreator
34 (2005,1137) PN [PDF_CONTROL_PREP_PARS] # 22, 1 Unknown Tag & Data
35 (2005,1138) PN (no value available) # 0, 0 Unknown Tag & Data
36 (2005,1139) PN [IEEE_PDF] # 8, 1 Unknown Tag & Data
37 (2005,1140) PN (no value available) # 0, 0 Unknown Tag & Data
38 (2005,1141) PN (no value available) # 0, 0 Unknown Tag & Data
39 (2005,1143) SL 7934 # 4, 1 Unknown Tag & Data
40 (2005,1144) OW 19b6\0000\005f\0000\1b2a\0000\00f3\0000\1eee\0000\0000\0000\0008... # 7934, 1 Unknown
Tag & Data
41 (2005,1147) CS [Y] # 2, 1 Unknown Tag & Data
42 (fffe,e00d) na (ItemDelimitationItem) # 0, 0 ItemDelimitationItem
43 ...
44 """
45
46 import sys
47 import gdcm
48
49 if __name__ == "__main__":
50
51     file1 = sys.argv[1]
52     file2 = sys.argv[2]
53
54     r = gdcm.Reader()
55     r.SetFileName( file1 )
56     if not r.Read():
57         sys.exit(1)
58
59     fg = gdcm.FileNameGenerator()
60     f = r.GetFile()
61     ds = f.GetDataSet()
62     tsis = gdcm.Tag(0x2005,0x1132) #
63     if ds.FindDataElement( tsis ):
64         sis = ds.GetDataElement( tsis )
65         #sqsis = sis.GetSequenceOfItems()
66         # GetValueAsSQ handle more cases
67         sqsis = sis.GetValueAsSQ()
68         if sqsis.GetNumberOfItems():
69             nitems = sqsis.GetNumberOfItems();
70             fg.SetNumberOfFileNames( nitems )
71             fg.SetPrefix( file2 )
72             if not fg.Generate():
73                 print "problem"
74                 sys.exit(1)
75             for i in range(0,nitems):
76                 item1 = sqsis.GetItem(i+1) # Item start at 1
77                 nestedds = item1.GetNestedDataSet()
78                 tprcs = gdcm.Tag(0x2005,0x1144) #
79                 if nestedds.FindDataElement( tprcs ):
80                     prcs = nestedds.GetDataElement( tprcs )
81                     bv = prcs.GetByteValue()

```

```
82         print bv
83         f = open( fg.GetFilename(i) , "w" )
84         f.write( bv.WriteBuffer() )
```

# Index

- ~ASN1
  - gdcmm::ASN1, [168](#)
- ~AnonymizeEvent
  - gdcmm::AnonymizeEvent, [153](#)
- ~Anonymizer
  - gdcmm::Anonymizer, [156](#)
- ~Attribute
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, [184](#)
- ~AudioCodec
  - gdcmm::AudioCodec, [194](#)
- ~BaseCompositeMessage
  - gdcmm::network::BaseCompositeMessage, [198](#)
- ~BasePDU
  - gdcmm::network::BasePDU, [200](#)
- ~BaseRootQuery
  - gdcmm::BaseRootQuery, [202](#)
- ~Bitmap
  - gdcmm::Bitmap, [212](#)
- ~BitmapToBitmapFilter
  - gdcmm::BitmapToBitmapFilter, [218](#)
- ~BoxRegion
  - gdcmm::BoxRegion, [220](#)
- ~ByteSwapFilter
  - gdcmm::ByteSwapFilter, [224](#)
- ~ByteValue
  - gdcmm::ByteValue, [227](#)
- ~CAPICryptographicMessageSyntax
  - gdcmm::CAPICryptographicMessageSyntax, [232](#)
- ~CSAHeader
  - gdcmm::CSAHeader, [270](#)
- ~Coder
  - gdcmm::Coder, [245](#)
- ~Command
  - gdcmm::Command, [250](#)
- ~CommandDataSet
  - gdcmm::CommandDataSet, [252](#)
- ~CryptoFactory
  - gdcmm::CryptoFactory, [261](#)
- ~CryptographicMessageSyntax
  - gdcmm::CryptographicMessageSyntax, [263](#)
- ~Curve
  - gdcmm::Curve, [281](#)
- ~DICOMDIRGenerator
  - gdcmm::DICOMDIRGenerator, [312](#)
- ~DataEvent
  - gdcmm::DataEvent, [294](#)
- ~DataSetEvent
  - gdcmm::DataSetEvent, [303](#)
- ~Decoder
  - gdcmm::Decoder, [305](#)
- ~Defs
  - gdcmm::Defs, [307](#)
- ~DeltaEncodingCodec
  - gdcmm::DeltaEncodingCodec, [310](#)
- ~DictConverter
  - gdcmm::DictConverter, [317](#)
- ~DictPrinter
  - gdcmm::DictPrinter, [322](#)
- ~Dicts
  - gdcmm::Dicts, [323](#)
- ~DirectionCosines
  - gdcmm::DirectionCosines, [327](#)
- ~Directory
  - gdcmm::Directory, [329](#)
- ~Dumper
  - gdcmm::Dumper, [334](#)
- ~Element
  - gdcmm::Element< TVR, VM::VM1\_n >, [340](#)
- ~Event
  - gdcmm::Event, [358](#)
- ~Exception
  - gdcmm::Exception, [360](#)
- ~FileAnonymizer
  - gdcmm::FileAnonymizer, [371](#)
- ~FileChangeTransferSyntax
  - gdcmm::FileChangeTransferSyntax, [374](#)
- ~FileDerivation
  - gdcmm::FileDerivation, [375](#)
- ~FileExplicitFilter
  - gdcmm::FileExplicitFilter, [378](#)
- ~FileNameEvent
  - gdcmm::FileNameEvent, [388](#)
- ~FileStreamer
  - gdcmm::FileStreamer, [394](#)
- ~FilenameGenerator
  - gdcmm::FilenameGenerator, [390](#)
- ~Global
  - gdcmm::Global, [405](#)
- ~GroupDict

- gdcmm::GroupDict, 407
- ~IconImageFilter
  - gdcmm::IconImageFilter, 409
- ~IconImageGenerator
  - gdcmm::IconImageGenerator, 411
- ~Image
  - gdcmm::Image, 415
- ~ImageApplyLookupTable
  - gdcmm::ImageApplyLookupTable, 419
- ~ImageChangePhotometricInterpretation
  - gdcmm::ImageChangePhotometricInterpretation, 421
- ~ImageChangePlanarConfiguration
  - gdcmm::ImageChangePlanarConfiguration, 424
- ~ImageChangeTransferSyntax
  - gdcmm::ImageChangeTransferSyntax, 427
- ~ImageCodec
  - gdcmm::ImageCodec, 431
- ~ImageConverter
  - gdcmm::ImageConverter, 436
- ~ImageFragmentSplitter
  - gdcmm::ImageFragmentSplitter, 439
- ~ImageReader
  - gdcmm::ImageReader, 445
- ~ImageRegionReader
  - gdcmm::ImageRegionReader, 448
- ~ImageToImageFilter
  - gdcmm::ImageToImageFilter, 451
- ~ImageWriter
  - gdcmm::ImageWriter, 453
- ~JPEG12Codec
  - gdcmm::JPEG12Codec, 473
- ~JPEG16Codec
  - gdcmm::JPEG16Codec, 475
- ~JPEG2000Codec
  - gdcmm::JPEG2000Codec, 477
- ~JPEG8Codec
  - gdcmm::JPEG8Codec, 481
- ~JPEGCodec
  - gdcmm::JPEGCodec, 484
- ~JPEGLSCodec
  - gdcmm::JPEGLSCodec, 488
- ~JSON
  - gdcmm::JSON, 491
- ~KAKADUCodec
  - gdcmm::KAKADUCodec, 493
- ~LookupTable
  - gdcmm::LookupTable, 498
- ~MD5
  - gdcmm::MD5, 505
- ~MemberCommand
  - gdcmm::MemberCommand, 515
- ~MeshPrimitive
  - gdcmm::MeshPrimitive, 519
- ~ModuleEntry
  - gdcmm::ModuleEntry, 524
- ~Object
  - gdcmm::Object, 536
- ~OpenSSLCryptographicMessageSyntax
  - gdcmm::OpenSSLCryptographicMessageSyntax, 539
- ~OpenSSL7CryptographicMessageSyntax
  - gdcmm::OpenSSL7CryptographicMessageSyntax, 543
- ~Orientation
  - gdcmm::Orientation, 545
- ~Overlay
  - gdcmm::Overlay, 549
- ~PDBHeader
  - gdcmm::PDBHeader, 562
- ~PDFCodec
  - gdcmm::PDFCodec, 564
- ~PGXCodec
  - gdcmm::PGXCodec, 568
- ~PNMCodec
  - gdcmm::PNMCodec, 590
- ~PVRGCodec
  - gdcmm::PVRGCodec, 611
- ~ParseException
  - gdcmm::ParseException, 553
- ~Parser
  - gdcmm::Parser, 555
- ~Pixmap
  - gdcmm::Pixmap, 578
- ~PixmapReader
  - gdcmm::PixmapReader, 582
- ~PixmapToPixmapFilter
  - gdcmm::PixmapToPixmapFilter, 584
- ~PixmapWriter
  - gdcmm::PixmapWriter, 587
- ~Preamble
  - gdcmm::Preamble, 592
- ~Printer
  - gdcmm::Printer, 603
- ~PrivateDict
  - gdcmm::PrivateDict, 605
- ~ProgressEvent
  - gdcmm::ProgressEvent, 609
- ~PythonFilter
  - gdcmm::PythonFilter, 613
- ~QueryBase
  - gdcmm::QueryBase, 614
- ~RAWCodec
  - gdcmm::RAWCodec, 626
- ~RLECodec
  - gdcmm::RLECodec, 638
- ~Reader
  - gdcmm::Reader, 630
- ~Region
  - gdcmm::Region, 633



- ~Rescaler
  - gdcm::Rescaler, [635](#)
- ~SHA1
  - gdcm::SHA1, [680](#)
- ~Scanner
  - gdcm::Scanner, [645](#)
- ~Segment
  - gdcm::Segment, [650](#)
- ~SegmentReader
  - gdcm::SegmentReader, [656](#)
- ~SegmentWriter
  - gdcm::SegmentWriter, [658](#)
- ~SegmentedPaletteColorLookupTable
  - gdcm::SegmentedPaletteColorLookupTable, [653](#)
- ~SerieHelper
  - gdcm::SerieHelper, [671](#)
- ~ServiceClassUser
  - gdcm::ServiceClassUser, [676](#)
- ~SimpleMemberCommand
  - gdcm::SimpleMemberCommand, [683](#)
- ~SimpleSubjectWatcher
  - gdcm::SimpleSubjectWatcher, [684](#)
- ~SmartPointer
  - gdcm::SmartPointer, [687](#)
- ~Sorter
  - gdcm::Sorter, [693](#)
- ~Spacing
  - gdcm::Spacing, [695](#)
- ~SplitMosaicFilter
  - gdcm::SplitMosaicFilter, [697](#)
- ~StreamImageReader
  - gdcm::StreamImageReader, [700](#)
- ~StreamImageWriter
  - gdcm::StreamImageWriter, [704](#)
- ~StringFilter
  - gdcm::StringFilter, [711](#)
- ~Subject
  - gdcm::Subject, [715](#)
- ~Surface
  - gdcm::Surface, [719](#)
- ~SurfaceReader
  - gdcm::SurfaceReader, [727](#)
- ~SurfaceWriter
  - gdcm::SurfaceWriter, [729](#)
- ~Table
  - gdcm::Table, [737](#)
- ~TableEntry
  - gdcm::TableEntry, [738](#)
- ~TableReader
  - gdcm::TableReader, [739](#)
- ~TableRow
  - gdcm::network::TableRow, [741](#)
- ~TagPath
  - gdcm::TagPath, [749](#)
- ~Testing
  - gdcm::Testing, [751](#)
- ~Trace
  - gdcm::Trace, [755](#)
- ~Transition
  - gdcm::network::Transition, [763](#)
- ~ULAction
  - gdcm::network::ULAction, [788](#)
- ~ULBasicCallback
  - gdcm::network::ULBasicCallback, [818](#)
- ~ULConnection
  - gdcm::network::ULConnection, [819](#)
- ~ULConnectionCallback
  - gdcm::network::ULConnectionCallback, [822](#)
- ~ULConnectionManager
  - gdcm::network::ULConnectionManager, [825](#)
- ~ULEvent
  - gdcm::network::ULEvent, [827](#)
- ~ULWritingCallback
  - gdcm::network::ULWritingCallback, [829](#)
- ~UserInformation
  - gdcm::network::UserInformation, [837](#)
- ~Validate
  - gdcm::Validate, [839](#)
- ~Value
  - gdcm::Value, [841](#)
- ~Version
  - gdcm::Version, [843](#)
- ~Writer
  - gdcm::Writer, [920](#)
- ~XMLDictReader
  - gdcm::XMLDictReader, [923](#)
- ~XMLPrinter
  - gdcm::XMLPrinter, [925](#)
- ~XMLPrivateDictReader
  - gdcm::XMLPrivateDictReader, [927](#)
- ~vtkGDCMImageReader
  - vtkGDCMImageReader, [861](#)
- ~vtkGDCMImageReader2
  - vtkGDCMImageReader2, [868](#)
- ~vtkGDCMImageWriter
  - vtkGDCMImageWriter, [873](#)
- ~vtkGDCMMedicalImageProperties
  - vtkGDCMMedicalImageProperties, [877](#)
- ~vtkGDCMPolyDataReader
  - vtkGDCMPolyDataReader, [879](#)
- ~vtkGDCMPolyDataWriter
  - vtkGDCMPolyDataWriter, [882](#)
- ~vtkGDCMTesting
  - vtkGDCMTesting, [885](#)
- ~vtkGDCMThreadedImageReader
  - vtkGDCMThreadedImageReader, [887](#)
- ~vtkGDCMThreadedImageReader2
  - vtkGDCMThreadedImageReader2, [890](#)

- ~vtkImageColorViewer
  - vtkImageColorViewer, [895](#)
- ~vtkImageMapToColors16
  - vtkImageMapToColors16, [901](#)
- ~vtkImageMapToWindowLevelColors2
  - vtkImageMapToWindowLevelColors2, [904](#)
- ~vtkImagePlanarComponentsToComponents
  - vtkImagePlanarComponentsToComponents, [906](#)
- ~vtkImageRGBToYBR
  - vtkImageRGBToYBR, [907](#)
- ~vtkImageYBRToRGB
  - vtkImageYBRToRGB, [909](#)
- ~vtkLookupTable16
  - vtkLookupTable16, [910](#)
- ~vtkRTStructSetProperties
  - vtkRTStructSetProperties, [913](#)
- AAAbortPDU
  - gdcm::network::AAAbortPDU, [140](#)
- AAAssociateACPDU
  - gdcm::network::AAAssociateACPDU, [143](#)
  - gdcm::network::AAAssociateRQPDU, [149](#)
- AAAssociateRJPDU
  - gdcm::network::AAAssociateRJPDU, [145](#)
- AAAssociateRQPDU
  - gdcm::network::AAAssociateACPDU, [143](#)
  - gdcm::network::AAAssociateRQPDU, [147](#)
- AE
  - gdcm::VR, [852](#)
- AEComp
  - gdcm, [124](#)
- AES128\_CIPHER
  - gdcm::CryptographicMessageSyntax, [262](#)
- AES192\_CIPHER
  - gdcm::CryptographicMessageSyntax, [262](#)
- AES256\_CIPHER
  - gdcm::CryptographicMessageSyntax, [263](#)
- ALGOType
  - gdcm::Segment, [650](#)
- ALGOType\_END
  - gdcm::Segment, [650](#)
- ARGB
  - gdcm::PhotometricInterpretation, [570](#)
- ARTIMTimer
  - gdcm::network::ARTIMTimer, [167](#)
- AReleaseRPPDU
  - gdcm::network::AReleaseRPPDU, [164](#)
- AReleaseRQPDU
  - gdcm::network::AReleaseRQPDU, [165](#)
- AS
  - gdcm::VR, [852](#)
- ASComp
  - gdcm, [124](#)
- ASN1
  - gdcm::ASN1, [168](#)
- AT
  - gdcm::VR, [852](#)
- AUTOMATIC
  - gdcm::Segment, [650](#)
- AXIAL
  - gdcm::Orientation, [545](#)
- AbstractSyntax
  - gdcm::network::AbstractSyntax, [151](#)
- ActiveComponent
  - vtkImageMapToColors16, [902](#)
- Add
  - gdcm::GroupDict, [407](#)
- AddAcceptedPresentationContext
  - gdcm::network::ULConnection, [819](#)
- AddCSAHeaderDictEntry
  - gdcm::CSAHeaderDict, [273](#)
- AddContourReferencedFrameOfReference
  - vtkRTStructSetProperties, [913](#)
- AddDerivationDescription
  - gdcm::FileDerivation, [375](#)
- AddDictEntry
  - gdcm::Dict, [314](#)
  - gdcm::PrivateDict, [605](#)
- AddFile
  - gdcm::FileSet, [391](#), [392](#)
  - gdcm::SerieHelper, [672](#)
- AddFileName
  - gdcm::SerieHelper, [672](#)
- AddFragment
  - gdcm::SequenceOfFragments, [662](#)
- AddGroupLength
  - gdcm::DictConverter, [317](#)
- AddIOD
  - gdcm::IODs, [463](#)
- AddIODEntry
  - gdcm::IOD, [460](#)
- AddImageDirectoryRecord
  - gdcm::DICOMDIRGenerator, [312](#)
- AddInput
  - vtkImageColorViewer, [895](#)
- AddInputConnection
  - vtkImageColorViewer, [895](#)
- AddItem
  - gdcm::SequenceOfItems, [667](#)
- AddMacro
  - gdcm::Macros, [503](#)
  - gdcm::Module, [522](#)
- AddMacroEntry
  - gdcm::Macro, [502](#)
- AddModule
  - gdcm::Modules, [526](#)
- AddModuleEntry
  - gdcm::Module, [522](#)

- gdcmm::NestedModuleEntries, 533
- AddObserver
  - gdcmm::Subject, 715
- AddPatientDirectoryRecord
  - gdcmm::DICOMDIRGenerator, 312
- AddPresentationContext
  - gdcmm::PresentationContextGenerator, 596
  - gdcmm::network::AAssociateRQPDU, 147
- AddPresentationContextAC
  - gdcmm::network::AAssociateACPDU, 143
- AddPresentationDataValue
  - gdcmm::network::PDataTFPDU, 558
- AddPrimitiveData
  - gdcmm::MeshPrimitive, 519
- AddPrivateTag
  - gdcmm::Scanner, 645
- AddPurposeOfReferenceCodeSequence
  - gdcmm::FileDerivation, 375
- AddQueryDataSet
  - gdcmm::BaseRootQuery, 202
- AddReference
  - gdcmm::FileDerivation, 375
- AddReferencedFrameOfReference
  - vtkRTStructSetProperties, 914
- AddRestriction
  - gdcmm::SerieHelper, 672
- AddRoleSelectionSub
  - gdcmm::network::UserInformation, 837
- AddSOPClassExtendedNegociationSub
  - gdcmm::network::UserInformation, 837
- AddSegment
  - gdcmm::SegmentWriter, 658
- AddSelect
  - gdcmm::Sorter, 693
- AddSeriesDirectoryRecord
  - gdcmm::DICOMDIRGenerator, 312
- AddSkipTag
  - gdcmm::Scanner, 645
- AddSourceImageSequence
  - gdcmm::FileDerivation, 376
- AddStructureSetROI
  - vtkRTStructSetProperties, 914
- AddStructureSetROIObservation
  - vtkRTStructSetProperties, 914
- AddStudyDirectoryRecord
  - gdcmm::DICOMDIRGenerator, 312
- AddSurface
  - gdcmm::Segment, 650
- AddTag
  - gdcmm::Scanner, 645
- AddTransferSyntax
  - gdcmm::PresentationContext, 593
  - gdcmm::network::PresentationContextRQ, 598
- AffectedSOPClassUID
  - gdcmm::network::CEchoRQ, 234
- Allocate
  - gdcmm::LookupTable, 498
- AmbulatoryECGWaveformStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 775
- AnatomicRegion
  - gdcmm::Segment, 651
- AnonymizeEvent
  - gdcmm::AnonymizeEvent, 153
- Anonymizer
  - gdcmm::Anonymizer, 156
- Append
  - gdcmm::ByteValue, 227
  - gdcmm::Global, 405
- AppendFrameEncode
  - gdcmm::ImageCodec, 431
  - gdcmm::JPEG2000Codec, 478
  - gdcmm::JPEGCodec, 484
  - gdcmm::JPEGLSCodec, 488
  - gdcmm::RLECodec, 638
- AppendImplementationClassUID
  - gdcmm::FileMetaInformation, 381
- AppendRowEncode
  - gdcmm::ImageCodec, 431
  - gdcmm::JPEG2000Codec, 478
  - gdcmm::JPEGCodec, 484
  - gdcmm::JPEGLSCodec, 488
  - gdcmm::RLECodec, 639
- AppendToDataElement
  - gdcmm::FileStreamer, 394
- AppendToGroupDataElement
  - gdcmm::FileStreamer, 394
- ApplicationContext
  - gdcmm::network::ApplicationContext, 160
- Apply
  - gdcmm::ImageApplyLookupTable, 419
- ApplyInverseVideo
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- ApplyLookupTable
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- ApplyPlanarConfiguration
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- ApplyShiftScale
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- ApplyYBRToRGB
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- AreOverlaysInPixelData
  - gdcmm::Bitmap, 212

- gdcm::Pixmap, [578](#)
- Area
  - gdcm::BoxRegion, [220](#)
  - gdcm::Region, [633](#)
- ArrayIncludeMacrosType
  - gdcm::Macro, [501](#)
  - gdcm::Module, [521](#)
- ArrayType
  - gdcm::Attribute, [171](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [177](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [184](#)
- AsynchronousOperationsWindowSub
  - gdcm::network::AsynchronousOperationsWindow↔ Sub, [168](#)
- Attribute
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [184](#)
  - gdcm::terminal, [137](#)
- Audio
  - gdcm::MediaStorage, [510](#)
- AudioCodec
  - gdcm::AudioCodec, [194](#)
- AudioSRStorageTrialRetired
  - gdcm::UIDs, [776](#)
- AutoPixelMinMax
  - gdcm::IconImageGenerator, [411](#)
- BALCPPProtect
  - gdcm::Anonymizer, [156](#)
- BLUE
  - gdcm::LookupTable, [498](#)
- BOOL\_FUNCTION\_PFILE\_PFILE\_POINTER
  - gdcm, [124](#)
- backslash
  - gdcm, [126](#)
- BadBigEndian
  - gdcm::SwapCode, [730](#)
- BadLittleEndian
  - gdcm::SwapCode, [730](#)
- BaseRootQuery
  - gdcm::BaseRootQuery, [202](#)
- BasicAnnotationBoxSOPClass
  - gdcm::UIDs, [774](#)
- BasicApplicationLevelConfidentialityProfile
  - gdcm::Anonymizer, [156](#)
- BasicCodedEntry
  - gdcm::SegmentHelper::BasicCodedEntry, [206](#)
- BasicColorImageBoxSOPClass
  - gdcm::UIDs, [774](#)
- BasicColorPrintManagementMetaSOPClass
  - gdcm::UIDs, [774](#)
- BasicFilmBoxSOPClass
  - gdcm::UIDs, [774](#)
- BasicFilmSessionSOPClass
  - gdcm::UIDs, [774](#)
- BasicGrayscaleImageBoxSOPClass
  - gdcm::UIDs, [774](#)
- BasicGrayscalePrintManagementMetaSOPClass
  - gdcm::UIDs, [774](#)
- BasicOffsetTable
  - gdcm::BasicOffsetTable, [208](#)
- BasicPrintImageOverlayBoxSOPClassRetired
  - gdcm::UIDs, [775](#)
- BasicStudyContentNotificationSOPClassRetired
  - gdcm::UIDs, [774](#)
- BasicTextSR
  - gdcm::MediaStorage, [509](#)
- BasicTextSRStorage
  - gdcm::UIDs, [776](#)
- BasicVoiceAudioWaveformStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [775](#)
- Begin
  - gdcm::CSAHeaderDict, [273](#)
  - gdcm::DataSet, [297](#)
  - gdcm::Dict, [314](#)
  - gdcm::IODs, [463](#)
  - gdcm::Scanner, [645](#)
  - gdcm::SequenceOfFragments, [662](#)
  - gdcm::SequenceOfItems, [667](#)
- BigEndian
  - gdcm::SwapCode, [730](#)
- BitSample
  - gdcm::JPEGCodec, [486](#)
  - gdcm::LookupTable, [500](#)
- Bitmap
  - gdcm::Bitmap, [212](#)
  - gdcm::JPEG2000Codec, [479](#)
  - gdcm::PixelFormat, [576](#)
- BitmapToBitmapFilter
  - gdcm::BitmapToBitmapFilter, [218](#)
- black
  - gdcm::terminal, [137](#)
- BlendingSoftcopyPresentationStateStorageSOPClass
  - gdcm::UIDs, [776](#)
- blink
  - gdcm::terminal, [137](#)
- blue
  - gdcm::terminal, [138](#)
- BoundingBox
  - gdcm::BoxRegion, [220](#)
- BoxRegion
  - gdcm::BoxRegion, [220](#)
- BreakConnection
  - gdcm::network::ULConnectionManager, [825](#)
- BreakConnectionNow

- gdcm::network::ULConnectionManager, 825
- BreastImagingRelevantPatientInformationQuery
  - gdcm::UIDs, 777
- BreastTomosynthesisImageStorage
  - gdcm::MediaStorage, 510
  - gdcm::UIDs, 779
- bright
  - gdcm::terminal, 137
- Build
  - vtkLookupTable16, 910
- ByteBuffer
  - gdcm::ByteBuffer, 222
- ByteSwap
  - gdcm::ByteSwapFilter, 224
- ByteSwapFilter
  - gdcm::ByteSwapFilter, 224
- ByteValue
  - gdcm::ByteValue, 226
- bytes
  - gdcm::Tag, 748
- C\_CANCEL\_RQ
  - gdcm::network::DIMSE, 326
- C\_ECHO\_RQ
  - gdcm::network::DIMSE, 326
- C\_ECHO\_RSP
  - gdcm::network::DIMSE, 326
- C\_FIND\_RQ
  - gdcm::network::DIMSE, 325
- C\_FIND\_RSP
  - gdcm::network::DIMSE, 325
- C\_GET\_RQ
  - gdcm::network::DIMSE, 325
- C\_GET\_RSP
  - gdcm::network::DIMSE, 325
- C\_MOVE\_RQ
  - gdcm::network::DIMSE, 325
- C\_MOVE\_RSP
  - gdcm::network::DIMSE, 326
- C\_STORE\_RQ
  - gdcm::network::DIMSE, 325
- C\_STORE\_RSP
  - gdcm::network::DIMSE, 325
- CALIBRATED
  - gdcm::Spacing, 695
- CAPI
  - gdcm::CryptoFactory, 261
- CAPICryptoFactory
  - gdcm::CAPICryptoFactory, 230
- CAPICryptographicMessageSyntax
  - gdcm::CAPICryptographicMessageSyntax, 232
- CEcho
  - gdcm::CompositeNetworkFunctions, 254
- CFind
  - gdcm::CompositeNetworkFunctions, 256
- CM
  - gdcm::SegmentHelper::BasicCodedEntry, 206
- CMYK
  - gdcm::PhotometricInterpretation, 570
- cMaxEventID
  - gdcm::network, 136
- cMaxStateID
  - gdcm::network, 136
- CMove
  - gdcm::CompositeNetworkFunctions, 256
- CONDENSED\_STYLE
  - gdcm::Printer, 603
- CONSOLE
  - gdcm::terminal, 138
- CORONAL
  - gdcm::Orientation, 545
- CS
  - gdcm::VR, 852
- CSAElement
  - gdcm::CSAElement, 266
- CSAHeader
  - gdcm::CSAHeader, 270
  - gdcm::DataSet, 301
- CSAHeaderDict
  - gdcm::CSAHeaderDict, 273
- CSAHeaderDictEntry
  - gdcm::CSAHeaderDictEntry, 275
- CSAHeaderType
  - gdcm::CSAHeader, 270
- CSANonImageStorage
  - gdcm::MediaStorage, 509
- CSComp
  - gdcm, 124
- CSD
  - gdcm::SegmentHelper::BasicCodedEntry, 206
- CSV
  - gdcm::SegmentHelper::BasicCodedEntry, 206
- CStore
  - gdcm::CompositeNetworkFunctions, 257
- CT\_private\_ELE
  - gdcm::TransferSyntax, 760
- CTImageStorage
  - gdcm::MediaStorage, 508
  - gdcm::UIDs, 775
- CV
  - gdcm::SegmentHelper::BasicCodedEntry, 206
- CanCode
  - gdcm::AudioCodec, 194
  - gdcm::Coder, 245
  - gdcm::ImageCodec, 431
  - gdcm::JPEG2000Codec, 478
  - gdcm::JPEGCodec, 484
  - gdcm::JPEGLSCCodec, 489

- gdcmm::KAKADUCodec, 493
- gdcmm::PDFCodec, 564
- gdcmm::PGXCodec, 568
- gdcmm::PNMCodec, 590
- gdcmm::PVRGCodec, 611
- gdcmm::RAWCodec, 626
- gdcmm::RLECodec, 639
- CanDecode
  - gdcmm::AudioCodec, 195
  - gdcmm::Decoder, 305
  - gdcmm::DeltaEncodingCodec, 310
  - gdcmm::ImageCodec, 432
  - gdcmm::JPEG2000Codec, 478
  - gdcmm::JPEGCodec, 484
  - gdcmm::JPEGLSCodec, 489
  - gdcmm::KAKADUCodec, 493
  - gdcmm::PDFCodec, 564
  - gdcmm::PGXCodec, 568
  - gdcmm::PNMCodec, 590
  - gdcmm::PVRGCodec, 611
  - gdcmm::RAWCodec, 626
  - gdcmm::RLECodec, 639
- CanDisplay
  - gdcmm::VR, 854
- CanEmptyTag
  - gdcmm::Anonymizer, 156
- CanRead
  - gdcmm::Reader, 630
- CanReadFile
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
- CanReadImage
  - gdcmm::StreamImageReader, 700
- CanStoreLossy
  - gdcmm::TransferSyntax, 760
- CanWriteFile
  - gdcmm::StreamImageWriter, 704
- CardiacElectrophysiologyWaveformStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 775
- CardiacRelevantPatientInformationQuery
  - gdcmm::UIDs, 778
- Change
  - gdcmm::FileChangeTransferSyntax, 374
  - gdcmm::FileExplicitFilter, 378
  - gdcmm::ImageChangePhotometricInterpretation, 421
  - gdcmm::ImageChangePlanarConfiguration, 424
  - gdcmm::ImageChangeTransferSyntax, 427
- ChangeFMI
  - gdcmm::FileExplicitFilter, 378
- ChangeMonochrome
  - gdcmm::ImageChangePhotometricInterpretation, 421
- CharacterDataHandler
  - gdcmm::TableReader, 739
- gdcmm::XMLDictReader, 923
- gdcmm::XMLPrivateDictReader, 927
- CheckDataElement
  - gdcmm::FileStreamer, 394
- CheckEvent
  - gdcmm::AnonymizeEvent, 153
  - gdcmm::DataEvent, 294
  - gdcmm::DataSetEvent, 303
  - gdcmm::Event, 358
  - gdcmm::FileNameEvent, 388
  - gdcmm::ProgressEvent, 609
- CheckFileMetaInformationOff
  - gdcmm::Writer, 920
- CheckFileMetaInformationOn
  - gdcmm::Writer, 920
- CheckTemplateFileName
  - gdcmm::FileStreamer, 394
- ChestCADSRStorage
  - gdcmm::UIDs, 777
- CipherTypes
  - gdcmm::CryptographicMessageSyntax, 262
- Clear
  - gdcmm::Bitmap, 212
  - gdcmm::ByteValue, 227
  - gdcmm::DataElement, 286
  - gdcmm::DataSet, 297
  - gdcmm::IOD, 460
  - gdcmm::IODs, 463
  - gdcmm::Item, 469
  - gdcmm::LookupTable, 498
  - gdcmm::Macro, 502
  - gdcmm::Macros, 503
  - gdcmm::Module, 522
  - gdcmm::Modules, 526
  - gdcmm::Preamble, 592
  - gdcmm::SequenceOfFragments, 662
  - gdcmm::SequenceOfItems, 667
  - gdcmm::SerieHelper, 672
  - gdcmm::Value, 841
  - vtkGDCMMedicalImageProperties, 877
  - vtkRTStructSetProperties, 914
- ClearInternalUIDs
  - gdcmm::Anonymizer, 157
- ClearSkipTags
  - gdcmm::Scanner, 646
- ClearTags
  - gdcmm::Scanner, 646
- Clone
  - gdcmm::BoxRegion, 221
  - gdcmm::ImageCodec, 432
  - gdcmm::JPEG2000Codec, 478
  - gdcmm::JPEGCodec, 484
  - gdcmm::JPEGLSCodec, 489
  - gdcmm::KAKADUCodec, 493



- gdcmm::PGXCodec, 568
- gdcmm::PNMCodec, 590
- gdcmm::PVRGCodec, 611
- gdcmm::RAWCodec, 626
- gdcmm::RLECodec, 639
- gdcmm::Region, 633
- Code
  - gdcmm::Coder, 245
  - gdcmm::JPEG2000Codec, 478
  - gdcmm::JPEGCodec, 484
  - gdcmm::JPEGLSCodec, 489
  - gdcmm::JSON, 491
  - gdcmm::KAKADUCodec, 493
  - gdcmm::PVRGCodec, 612
  - gdcmm::RAWCodec, 626
  - gdcmm::RLECodec, 639
- CodeString
  - gdcmm::CodeString, 247
- Color
  - gdcmm::terminal, 137
- ColorArray
  - gdcmm::SurfaceHelper, 723
- ColorSoftcopyPresentationStateStorageSOPClass
  - gdcmm::UIDs, 776
- Command
  - gdcmm::Command, 250
- CommandDataSet
  - gdcmm::CommandDataSet, 252
- CommandTypes
  - gdcmm::network::DIMSE, 325
- CompOperators
  - gdcmm, 124
- Compatible
  - gdcmm::VM, 849
  - gdcmm::VR, 854
- Component
  - gdcmm::PersonName, 566
- ComprehensiveSR
  - gdcmm::MediaStorage, 509
- ComprehensiveSRStorage
  - gdcmm::UIDs, 776
- ComprehensiveSRStorageTrialRetired
  - gdcmm::UIDs, 776
- CompressionTypes
  - vtkGDCMImageWriter, 873
- Compute
  - gdcmm::MD5, 505
  - gdcmm::SHA1, 680
- ComputeBoundingBox
  - gdcmm::BoxRegion, 221
  - gdcmm::Region, 633
- ComputeBufferLength
  - gdcmm::ImageRegionReader, 448
- ComputeByteLength
  - gdcmm::SequenceOfFragments, 662
- ComputeDataElement
  - gdcmm::DataSet, 297
- ComputeDataSetMediaStorageSOPClass
  - gdcmm::FileMetaInformation, 381
- ComputeDataSetTransferSyntax
  - gdcmm::FileMetaInformation, 381
- ComputeDistAlongNormal
  - gdcmm::DirectionCosines, 327
- ComputeFile
  - gdcmm::MD5, 505
  - gdcmm::SHA1, 680
- ComputeFileMD5
  - gdcmm::Testing, 751
- ComputeGroupLength
  - gdcmm::DataSet, 298
- ComputeInterceptSlopePixelType
  - gdcmm::Rescaler, 635
- ComputeLength
  - gdcmm::ByteValue, 227
  - gdcmm::Fragment, 403
  - gdcmm::SequenceOfFragments, 662
  - gdcmm::SequenceOfItems, 667
- ComputeLossyFlag
  - gdcmm::Bitmap, 212
- ComputeMD5
  - gdcmm::Testing, 751
- ComputeMOSAICDimensions
  - gdcmm::SplitMosaicFilter, 697
- ComputeMediaStorageFromModality
  - gdcmm::ImageHelper, 440
- ComputeNumberOfSurfaces
  - gdcmm::SurfaceWriter, 729
- ComputeOffsetTable
  - gdcmm::JPEGCodec, 485
- ComputePixelAspectRatioFromPixelSpacing
  - gdcmm::Spacing, 695
- ComputePixelTypeFromMinMax
  - gdcmm::Rescaler, 635
- ComputeSpacingFromImagePositionPatient
  - gdcmm::ImageHelper, 440
- ComputeVR
  - gdcmm::DataSetHelper, 304
- ComputeZSpacing
  - gdcmm::IPPSorter, 467
- ComputedRadiographyImageStorage
  - gdcmm::MediaStorage, 508
  - gdcmm::UIDs, 775
- ConcatenatePDVBlobs
  - gdcmm::network::PresentationDataValue, 600
- ConcatenatePDVBlobsAsExplicit
  - gdcmm::network::PresentationDataValue, 600
- Conditional
  - gdcmm::Usage, 835

- const
  - gdcm::SOPClassUIDToIOD, 690
- const\_iterator
  - gdcm::CodeString, 247
  - gdcm::LO, 495
  - gdcm::String, 709
- const\_reference
  - gdcm::CodeString, 247
  - gdcm::LO, 495
  - gdcm::String, 709
- const\_reverse\_iterator
  - gdcm::CodeString, 247
  - gdcm::LO, 495
  - gdcm::String, 709
- ConstCharWrapper
  - gdcm::ConstCharWrapper, 258
- ConstIterator
  - gdcm::CSAHeaderDict, 273
  - gdcm::DataSet, 297
  - gdcm::Dict, 314
  - gdcm::Scanner, 645
  - gdcm::SequenceOfFragments, 661
  - gdcm::SequenceOfItems, 667
- Construct
  - gdcm::BaseRootQuery, 203
- ConstructAbortPDU
  - gdcm::network::PDUFactory, 565
- ConstructCEchoRQ
  - gdcm::network::CompositeMessageFactory, 253
- ConstructCFindRQ
  - gdcm::network::CompositeMessageFactory, 253
- ConstructCMoveRQ
  - gdcm::network::CompositeMessageFactory, 253
- ConstructCStoreRQ
  - gdcm::network::CompositeMessageFactory, 253
- ConstructCStoreRSP
  - gdcm::network::CompositeMessageFactory, 253
- ConstructFromString
  - gdcm::TagPath, 749
- ConstructFromTagList
  - gdcm::TagPath, 749
- ConstructPDU
  - gdcm::network::PDUFactory, 565
- ConstructPDV
  - gdcm::network::BaseCompositeMessage, 198
  - gdcm::network::CEchoRQ, 234
  - gdcm::network::CFindRQ, 238
  - gdcm::network::CMoveRQ, 241
  - gdcm::network::CStoreRQ, 278
  - gdcm::network::CStoreRSP, 279
- ConstructPDVByDataSet
  - gdcm::network::CEchoRSP, 235
  - gdcm::network::CFindCancelRQ, 236
  - gdcm::network::CFindRSP, 239
  - gdcm::network::CMoveCancelRq, 240
  - gdcm::network::CMoveRSP, 243
- ConstructQuery
  - gdcm::CompositeNetworkFunctions, 257
- ConstructReleasePDU
  - gdcm::network::PDUFactory, 565
- ConstructorType
  - gdcm::Dicts, 323
- Convert
  - gdcm::DictConverter, 317
  - gdcm::ImageConverter, 436
- ConvertRGBToPaletteColor
  - gdcm::IconImageGenerator, 411
- ConvertToCXX
  - gdcm::DictConverter, 317
- ConvertToXML
  - gdcm::DictConverter, 317
- Create
  - gdcm::Preamble, 592
- CreateCEchoPDU
  - gdcm::network::PDUFactory, 565
- CreateCFindPDU
  - gdcm::network::PDUFactory, 565
- CreateCMSProvider
  - gdcm::CAPICryptoFactory, 230
  - gdcm::CryptoFactory, 261
  - gdcm::OpenSSLCryptoFactory, 538
  - gdcm::OpenSSLP7CryptoFactory, 541
- CreateCMovePDU
  - gdcm::network::PDUFactory, 565
- CreateCStoreRQPDU
  - gdcm::network::PDUFactory, 565
- CreateCStoreRSPPDU
  - gdcm::network::PDUFactory, 565
- CreateDefaultUniqueSeriesIdentifier
  - gdcm::SerieHelper, 672
- CreateUniqueSeriesIdentifier
  - gdcm::SerieHelper, 672
- Cross
  - gdcm::DirectionCosines, 327
- CrossDot
  - gdcm::DirectionCosines, 327
- CryptoFactory
  - gdcm::CryptoFactory, 261
- CryptoLib
  - gdcm::CryptoFactory, 261
- CryptographicMessageSyntax
  - gdcm::CryptographicMessageSyntax, 263
- Curve
  - gdcm::Curve, 281
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
- Curves
  - gdcm::Pixmap, 579



- cyan
  - gdcm::terminal, [138](#)
- DA
  - gdcm::VR, [852](#)
- DAComp
  - gdcm, [124](#)
- DATASET\_FORMAT
  - gdcm::CSAHeader, [270](#)
- DEFAULT
  - gdcm::CryptoFactory, [261](#)
- DES3\_CIPHER
  - gdcm::CryptographicMessageSyntax, [262](#)
- DETECTOR
  - gdcm::Spacing, [695](#)
- DICOMApplicationContextName
  - gdcm::UIDs, [774](#)
- DICOMControlledTerminology
  - gdcm::UIDs, [774](#)
- DICOMDIR
  - gdcm::DICOMDIR, [311](#)
- DICOMDIRGenerator
  - gdcm::DICOMDIRGenerator, [312](#)
- DICOMUIDRegistry
  - gdcm::UIDs, [774](#)
- DICT\_DEBUG
  - gdcm::DictConverter, [317](#)
- DICT\_DEFAULT
  - gdcm::DictConverter, [317](#)
- DICT\_XML
  - gdcm::DictConverter, [317](#)
- DS
  - gdcm::VR, [853](#)
- DT
  - gdcm::VR, [853](#)
- DTComp
  - gdcm, [124](#)
- DataElement
  - gdcm::DataElement, [286](#)
  - gdcm::Value, [842](#)
- DataElementSet
  - gdcm::DataSet, [297](#)
- DataElementType
  - gdcm::ModuleEntry, [525](#)
- DataEvent
  - gdcm::DataEvent, [293](#), [294](#)
- DataField
  - gdcm::CSAElement, [268](#)
- DataPtr
  - gdcm::CSAElement, [266](#)
- DataSetEvent
  - gdcm::DataSetEvent, [303](#)
- DataSetHandled
  - gdcm::network::ULConnectionCallback, [822](#)
- DataSetHandles
  - gdcm::network::ULConnectionCallback, [822](#)
- DataSetMS
  - gdcm::FileMetaInformation, [383](#)
- DataSetTS
  - gdcm::FileMetaInformation, [383](#)
- DataWasPassed
  - vtkImageMapToColors16, [902](#)
- DebugOff
  - gdcm::Trace, [755](#)
- DebugOn
  - gdcm::Trace, [755](#)
- Decode
  - gdcm::AudioCodec, [195](#)
  - gdcm::Base64, [195](#)
  - gdcm::Curve, [281](#)
  - gdcm::Decoder, [305](#)
  - gdcm::DeltaEncodingCodec, [310](#)
  - gdcm::ImageCodec, [432](#)
  - gdcm::JPEG2000Codec, [478](#)
  - gdcm::JPEGCodec, [485](#)
  - gdcm::JPEGLSCodec, [489](#)
  - gdcm::JSON, [491](#)
  - gdcm::KAKADUCodec, [493](#)
  - gdcm::LookupTable, [498](#)
  - gdcm::PDFCodec, [564](#)
  - gdcm::PVRGCodec, [612](#)
  - gdcm::RAWCodec, [626](#)
  - gdcm::RLECodec, [639](#)
- DecodeByStreams
  - gdcm::Decoder, [305](#)
  - gdcm::ImageCodec, [432](#)
  - gdcm::JPEG12Codec, [473](#)
  - gdcm::JPEG16Codec, [475](#)
  - gdcm::JPEG2000Codec, [478](#)
  - gdcm::JPEG8Codec, [481](#)
  - gdcm::JPEGCodec, [485](#)
  - gdcm::RAWCodec, [626](#)
  - gdcm::RLECodec, [639](#)
- DecodeBytes
  - gdcm::RAWCodec, [627](#)
- DecodeExtent
  - gdcm::JPEG2000Codec, [478](#)
  - gdcm::JPEGCodec, [485](#)
  - gdcm::JPEGLSCodec, [489](#)
  - gdcm::RLECodec, [639](#)
- Decompress
  - gdcm::Overlay, [549](#)
- Decrypt
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [539](#)
  - gdcm::OpenSSLP7CryptographicMessageSyntax, [543](#)

- DeepCopy
  - vtkRTStructSetProperties, [914](#)
- Default
  - gdcm::FileMetaInformation, [381](#)
- DefinePixelExtent
  - gdcm::StreamImageReader, [700](#)
  - gdcm::StreamImageWriter, [705](#)
- DefineProperBufferLength
  - gdcm::StreamImageReader, [701](#)
  - gdcm::StreamImageWriter, [705](#)
- DefinedTerms
  - gdcm::DefinedTerms, [306](#)
- DeflatedExplicitVRLittleEndian
  - gdcm::TransferSyntax, [759](#)
  - gdcm::UIDs, [772](#)
- DeformableSpatialRegistrationStorage
  - gdcm::UIDs, [776](#)
- Defs
  - gdcm::Defs, [307](#)
- DeleteDirectory
  - gdcm::System, [733](#)
- DeltaEncodingCodec
  - gdcm::DeltaEncodingCodec, [310](#)
- Derive
  - gdcm::FileDerivation, [376](#)
- Description
  - gdcm::ModuleEntry, [524](#)
- DescriptionField
  - gdcm::ModuleEntry, [525](#)
- DetachedInterpretationManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedPatientManagementMetaSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedPatientManagementSOPClass
  - gdcm::MediaStorage, [509](#)
- DetachedPatientManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedResultsManagementMetaSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedResultsManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedStudyManagementMetaSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedStudyManagementSOPClass
  - gdcm::MediaStorage, [509](#)
- DetachedStudyManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetachedVisitManagementSOPClass
  - gdcm::MediaStorage, [509](#)
- DetachedVisitManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- DetailSRStorageTrialRetired
  - gdcm::UIDs, [776](#)
- DetermineEventByPDU
  - gdcm::network::PDUFactory, [565](#)
- dicomAETitle
  - gdcm::UIDs, [778](#)
- dicomApplicationCluster
  - gdcm::UIDs, [778](#)
- dicomAssociationAcceptor
  - gdcm::UIDs, [778](#)
- dicomAssociationInitiator
  - gdcm::UIDs, [778](#)
- dicomAuthorizedNodeCertificateReference
  - gdcm::UIDs, [778](#)
- dicomConfigurationRoot
  - gdcm::UIDs, [778](#)
- dicomDescription
  - gdcm::UIDs, [778](#)
- dicomDevice
  - gdcm::UIDs, [778](#)
- dicomDeviceName
  - gdcm::UIDs, [778](#)
- dicomDeviceSerialNumber
  - gdcm::UIDs, [778](#)
- dicomDevicesRoot
  - gdcm::UIDs, [778](#)
- dicomHostname
  - gdcm::UIDs, [778](#)
- dicomInstalled
  - gdcm::UIDs, [778](#)
- dicomInstitutionAddress
  - gdcm::UIDs, [778](#)
- dicomInstitutionDepartmentName
  - gdcm::UIDs, [778](#)
- dicomInstitutionName
  - gdcm::UIDs, [778](#)
- dicomIssuerOfPatientID
  - gdcm::UIDs, [778](#)
- dicomManufacturer
  - gdcm::UIDs, [778](#)
- dicomManufacturerModelName
  - gdcm::UIDs, [778](#)
- dicomNetworkAE
  - gdcm::UIDs, [778](#)
- dicomNetworkConnection
  - gdcm::UIDs, [779](#)
- dicomNetworkConnectionReference
  - gdcm::UIDs, [778](#)
- dicomPort
  - gdcm::UIDs, [778](#)
- dicomPreferredCalledAETitle
  - gdcm::UIDs, [778](#)
- dicomPreferredCallingAETitle
  - gdcm::UIDs, [778](#)
- dicomPrimaryDeviceType
  - gdcm::UIDs, [778](#)
- dicomRelatedDeviceReference

- gdcmm::UIDs, [778](#)
- dicomSOPClass
  - gdcmm::UIDs, [778](#)
- dicomSoftwareVersion
  - gdcmm::UIDs, [778](#)
- dicomStationName
  - gdcmm::UIDs, [778](#)
- dicomSupportedCharacterSet
  - gdcmm::UIDs, [778](#)
- dicomTLSCyphersuite
  - gdcmm::UIDs, [778](#)
- dicomThisNodeCertificateReference
  - gdcmm::UIDs, [778](#)
- dicomTransferCapability
  - gdcmm::UIDs, [779](#)
- dicomTransferRole
  - gdcmm::UIDs, [778](#)
- dicomTransferSyntax
  - gdcmm::UIDs, [778](#)
- dicomUniqueAETitle
  - gdcmm::UIDs, [779](#)
- dicomUniqueAETitlesRegistryRoot
  - gdcmm::UIDs, [778](#)
- dicomVendorData
  - gdcmm::UIDs, [778](#)
- Dict
  - gdcmm::Dict, [314](#)
  - gdcmm::DictEntry, [320](#)
- DictConverter
  - gdcmm::DictConverter, [317](#)
- DictEntry
  - gdcmm::DictEntry, [319](#)
- DictPrinter
  - gdcmm::DictPrinter, [322](#)
- Dicts
  - gdcmm::CSAHeaderDict, [274](#)
  - gdcmm::Dict, [315](#)
  - gdcmm::Dicts, [323](#)
  - gdcmm::PrivateDict, [605](#)
- difference\_type
  - gdcmm::CodeString, [247](#)
  - gdcmm::LO, [495](#)
  - gdcmm::String, [709](#)
- DigitalIntraoralXRayImageStorageForPresentation
  - gdcmm::UIDs, [775](#)
- DigitalIntraoralXRayImageStorageForProcessing
  - gdcmm::MediaStorage, [508](#)
  - gdcmm::UIDs, [775](#)
- DigitalIntraoralXrayImageStorageForPresentation
  - gdcmm::MediaStorage, [508](#)
- DigitalMammographyImageStorageForPresentation
  - gdcmm::MediaStorage, [508](#)
- DigitalMammographyImageStorageForProcessing
  - gdcmm::MediaStorage, [508](#)
- DigitalMammographyXRayImageStorageForPresentation
  - gdcmm::UIDs, [775](#)
- DigitalMammographyXRayImageStorageForProcessing
  - gdcmm::UIDs, [775](#)
- DigitalXRayImageStorageForPresentation
  - gdcmm::MediaStorage, [508](#)
  - gdcmm::UIDs, [775](#)
- DigitalXRayImageStorageForProcessing
  - gdcmm::MediaStorage, [508](#)
  - gdcmm::UIDs, [775](#)
- dim
  - gdcmm::terminal, [137](#)
- Dimensions
  - gdcmm::Bitmap, [216](#)
  - gdcmm::ImageCodec, [435](#)
- DirCosTolerance
  - gdcmm::IPPSorter, [467](#)
- DirectionCosines
  - gdcmm::DirectionCosines, [327](#)
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- Directory
  - gdcmm::Directory, [329](#)
- DoByteSwap
  - gdcmm::ImageCodec, [432](#)
- DolconImage
  - gdcmm::PixmapWriter, [587](#)
- DoInvertMonochrome
  - gdcmm::ImageCodec, [432](#)
- DoOverlayCleanup
  - gdcmm::ImageCodec, [432](#)
- DoPaddedCompositePixelCode
  - gdcmm::ImageCodec, [432](#)
- DoPlanarConfiguration
  - gdcmm::ImageCodec, [432](#)
- DoSimpleCopy
  - gdcmm::ImageCodec, [432](#)
- DoYBR
  - gdcmm::ImageCodec, [432](#)
- Dot
  - gdcmm::DirectionCosines, [327](#)
- DropDuplicatePositions
  - gdcmm::IPPSorter, [467](#)
- Dumper
  - gdcmm::Dumper, [334](#)
- DuplicateAttributeError
  - gdcmm::Parser, [555](#)
- eAABORTPDUReturnedOpen
  - gdcmm::network, [135](#)
- eAABORTRequest
  - gdcmm::network, [135](#)
- eAASSOCIATE\_RQPDUReturned
  - gdcmm::network, [135](#)

- eAASSOCIATERequestLocalUser
  - gdcm::network, [135](#)
- eAASSOCIATEResponseAccept
  - gdcm::network, [135](#)
- eAASSOCIATEResponseReject
  - gdcm::network, [135](#)
- eARELEASE\_RPPDUReceived
  - gdcm::network, [135](#)
- eARELEASE\_RQPDUReceivedOpen
  - gdcm::network, [135](#)
- eARELEASERequest
  - gdcm::network, [135](#)
- eARELEASEResponse
  - gdcm::network, [135](#)
- eARTIMTimerExpired
  - gdcm::network, [135](#)
- eASSOCIATE\_ACPDUreceived
  - gdcm::network, [135](#)
- eASSOCIATE\_RJPDUreceived
  - gdcm::network, [135](#)
- eArabic
  - gdcm, [125](#)
- ECharSet
  - gdcm, [125](#)
- eCyrillic
  - gdcm, [125](#)
- EDGE
  - gdcm::MeshPrimitive, [518](#)
- eEventDoesNotExist
  - gdcm::network, [135](#)
- EEventID
  - gdcm::network, [135](#)
- eFind
  - gdcm, [125](#)
- eGB18030
  - gdcm, [125](#)
- eGreek
  - gdcm, [125](#)
- eHebrew
  - gdcm, [125](#)
- elImage
  - gdcm, [125](#)
- eJapanese
  - gdcm, [125](#)
- eJapaneseKanjiMultibyte
  - gdcm, [125](#)
- eJapaneseSupplementaryKanjiMultibyte
  - gdcm, [125](#)
- eKoreanHangulHanjaMultibyte
  - gdcm, [125](#)
- eLatin1
  - gdcm, [125](#)
- eLatin2
  - gdcm, [125](#)
- eLatin3
  - gdcm, [125](#)
- eLatin4
  - gdcm, [125](#)
- eLatin5
  - gdcm, [125](#)
- eMove
  - gdcm, [125](#)
- ePDATATFPDU
  - gdcm::network, [135](#)
- ePDATArequest
  - gdcm::network, [135](#)
- ePatient
  - gdcm, [125](#)
- ePatientRootType
  - gdcm, [126](#)
- EQueryLevel
  - gdcm, [125](#)
- EQueryType
  - gdcm, [125](#)
- ERootType
  - gdcm, [125](#)
- eSeries
  - gdcm, [125](#)
- eSta10ReleaseCollisionAc
  - gdcm::network, [136](#)
- eSta11ReleaseCollisionRq
  - gdcm::network, [136](#)
- eSta12ReleaseCollisionAcLocal
  - gdcm::network, [136](#)
- eSta13AwaitingClose
  - gdcm::network, [136](#)
- eSta1Idle
  - gdcm::network, [135](#)
- eSta2Open
  - gdcm::network, [135](#)
- eSta3WaitLocalAssoc
  - gdcm::network, [135](#)
- eSta4LocalAssocDone
  - gdcm::network, [135](#)
- eSta5WaitRemoteAssoc
  - gdcm::network, [135](#)
- eSta6TransferReady
  - gdcm::network, [135](#)
- eSta7WaitRelease
  - gdcm::network, [135](#)
- eSta8WaitLocalRelease
  - gdcm::network, [136](#)
- eSta9ReleaseCollisionRqLocal
  - gdcm::network, [136](#)
- eStaDoesNotExist
  - gdcm::network, [135](#)
- EStateID
  - gdcm::network, [135](#)

- eStudy
  - gdcm, [125](#)
- eStudyRootType
  - gdcm, [126](#)
- eThai
  - gdcm, [125](#)
- eTransportConnConfirmLocal
  - gdcm::network, [135](#)
- eTransportConnIndicLocal
  - gdcm::network, [135](#)
- eTransportConnectionClosed
  - gdcm::network, [135](#)
- eUTF8
  - gdcm, [125](#)
- eUnrecognizedPDURceived
  - gdcm::network, [135](#)
- elem
  - gdcm::SerieHelper::Rule, [642](#)
- Element
  - gdcm::Element< TVR, VM::VM1\_n >, [340](#)
- Empty
  - gdcm::Anonymizer, [157](#)
  - gdcm::BoxRegion, [221](#)
  - gdcm::DataElement, [286](#)
  - gdcm::FileAnonymizer, [371](#)
  - gdcm::Region, [633](#)
- EncapsulatedCDASStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [777](#)
- EncapsulatedDocument
  - gdcm::EncapsulatedDocument, [352](#)
- EncapsulatedPDFStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [777](#)
- Encode
  - gdcm::Base64, [196](#)
- EncodeBuffer
  - gdcm::JPEG12Codec, [473](#)
  - gdcm::JPEG16Codec, [475](#)
  - gdcm::JPEG8Codec, [481](#)
  - gdcm::JPEGCodec, [485](#)
- EncodeBytes
  - gdcm::System, [733](#)
- Encrypt
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [539](#)
  - gdcm::OpenSSL7CryptographicMessageSyntax, [543](#)
- End
  - gdcm::CSAHeaderDict, [274](#)
  - gdcm::DataSet, [298](#)
  - gdcm::Dict, [314](#)
  - gdcm::IODs, [463](#)
  - gdcm::Scanner, [646](#)
  - gdcm::SequenceOfFragments, [662](#)
  - gdcm::SequenceOfItems, [667](#)
- EndElement
  - gdcm::TableReader, [739](#)
  - gdcm::XMLDictReader, [923](#)
  - gdcm::XMLPrivateDictReader, [927](#)
- EndElementHandler
  - gdcm::Parser, [555](#)
- EndFilter
  - gdcm::SimpleSubjectWatcher, [684](#)
- EndWith
  - gdcm::Filename, [385](#)
- EnhancedCTImageStorage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- EnhancedMRImageStorage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- EnhancedSR
  - gdcm::MediaStorage, [509](#)
- EnhancedSRStorage
  - gdcm::UIDs, [776](#)
- EnhancedUSVolumeStorage
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [779](#)
- EnhancedXAImageStorage
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [776](#)
- EnhancedXRFIImageStorage
  - gdcm::UIDs, [776](#)
- EnumeratedValues
  - gdcm::EnumeratedValues, [357](#)
- ErrorOff
  - gdcm::Trace, [755](#)
- ErrorOn
  - gdcm::Trace, [755](#)
- ErrorType
  - gdcm::Parser, [555](#)
- EstablishConnection
  - gdcm::network::ULConnectionManager, [825](#)
- EstablishConnectionMove
  - gdcm::network::ULConnectionManager, [825](#)
- Event
  - gdcm::Event, [358](#)
- Exception
  - gdcm::Exception, [360](#)
- Execute
  - gdcm::Command, [250](#)
  - gdcm::MemberCommand, [515](#)
  - gdcm::SimpleMemberCommand, [683](#)
- ExecuteData
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMThreadedImageReader, [888](#)

- ExecuteInformation
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMThreadedImageReader, [888](#)
- ExecuteQuery
  - gdcm::StringFilter, [711](#)
- Explicit
  - gdcm::TransferSyntax, [759](#)
- ExplicitVRBigEndian
  - gdcm::TransferSyntax, [759](#)
  - gdcm::UIDs, [772](#)
- ExplicitVRLittleEndian
  - gdcm::TransferSyntax, [759](#)
  - gdcm::UIDs, [772](#)
- Explore
  - gdcm::Directory, [329](#)
- Extract
  - gdcm::IconImageFilter, [409](#)
- ExtractIconImages
  - gdcm::IconImageFilter, [409](#)
- ExtractVeprolIconImages
  - gdcm::IconImageFilter, [409](#)
- F
  - gdcm::Printer, [604](#)
  - gdcm::Reader, [632](#)
  - gdcm::Validate, [840](#)
  - gdcm::XMLPrinter, [925](#)
- FACET
  - gdcm::MeshPrimitive, [518](#)
- FD
  - gdcm::VR, [853](#)
- FL
  - gdcm::VR, [853](#)
- FLOAT16
  - gdcm::PixelFormat, [574](#)
- FLOAT32
  - gdcm::PixelFormat, [574](#)
- FLOAT64
  - gdcm::PixelFormat, [574](#)
- Fiducials
  - gdcm::Fiducials, [365](#)
- File
  - gdcm::File, [367](#)
- FileAnonymizer
  - gdcm::FileAnonymizer, [371](#)
- FileChangeTransferSyntax
  - gdcm::FileChangeTransferSyntax, [374](#)
  - gdcm::ImageCodec, [434](#)
- FileDerivation
  - gdcm::FileDerivation, [375](#)
- FileExists
  - gdcm::System, [733](#)
- FileExplicitFilter
  - gdcm::FileExplicitFilter, [378](#)
- FilesDirectory
  - gdcm::System, [734](#)
- FilesSymlink
  - gdcm::System, [734](#)
- FileList
  - gdcm, [124](#)
- FileMetaInformation
  - gdcm::FileMetaInformation, [381](#)
- FileName
  - vtkGDCMPolyDataReader, [880](#)
- FileNameEvent
  - gdcm::FileNameEvent, [388](#)
- FileNameOrdering
  - gdcm::SerieHelper, [672](#)
- FileNames
  - vtkGDCMImageReader, [864](#)
- FileSet
  - gdcm::FileSet, [391](#)
- FileSize
  - gdcm::System, [734](#)
- FileStreamer
  - gdcm::FileStreamer, [394](#)
- FileTime
  - gdcm::System, [734](#)
- FileType
  - gdcm::FileSet, [391](#)
- FileWithName
  - gdcm::FileWithName, [397](#)
- Filename
  - gdcm::Filename, [385](#)
- filename
  - gdcm::FileWithName, [397](#)
- FilenameGenerator
  - gdcm::FilenameGenerator, [389](#)
- FilenameType
  - gdcm::DICOMDIRGenerator, [312](#)
  - gdcm::Directory, [329](#)
  - gdcm::FilenameGenerator, [389](#)
- Filenames
  - gdcm::Sorter, [694](#)
- FilenamesType
  - gdcm::DICOMDIRGenerator, [312](#)
  - gdcm::Directory, [329](#)
  - gdcm::FilenameGenerator, [389](#)
- FilesType
  - gdcm::FileSet, [391](#)
- Fill
  - gdcm::ByteValue, [227](#)
- FillFromDataSet
  - gdcm::FileMetaInformation, [381](#)
- FillMedicalImageInformation
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [868](#)
  - vtkGDCMPolyDataReader, [879](#)

- FindCSAElementByName
  - gdcm::CSAHeader, [270](#)
- FindContext
  - gdcm::network::ULConnection, [819](#)
- FindDataElement
  - gdcm::DataSet, [298](#)
  - gdcm::Item, [469](#)
  - gdcm::SequenceOfItems, [668](#)
- FindDictEntry
  - gdcm::PrivateDict, [605](#)
- FindMacroEntry
  - gdcm::Macro, [502](#)
- FindModuleEntryInMacros
  - gdcm::Module, [522](#)
- FindNextDataElement
  - gdcm::DataSet, [298](#)
- FindPDBelementByName
  - gdcm::PDBHeader, [562](#)
- FindPatientRootQuery
  - gdcm::FindPatientRootQuery, [398](#)
- FindStudyRootQuery
  - gdcm::FindStudyRootQuery, [401](#)
- FirstRender
  - vtkImageColorViewer, [898](#)
- ForceRescale
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- FormatDateTime
  - gdcm::System, [734](#)
- Fragment
  - gdcm::Fragment, [403](#)
- FragmentVector
  - gdcm::SequenceOfFragments, [661](#)
- FromString
  - gdcm::StringFilter, [712](#)
- FujiPrivateCRImageStorage
  - gdcm::MediaStorage, [510](#)
- GDCM\_DIFFERENT
  - gdcm, [124](#)
- GDCM\_DO\_JOIN
  - gdcmStaticAssert.h, [1134](#)
- GDCM\_DO\_JOIN2
  - gdcmStaticAssert.h, [1134](#)
- GDCM\_EQUAL
  - gdcm, [124](#)
- GDCM\_EXPORT
  - gdcmWin32.h, [1188](#)
- GDCM\_FUNCTION
  - gdcmTrace.h, [1155](#)
- GDCM\_GREATER
  - gdcm, [125](#)
- GDCM\_GREATEROREQUAL
  - gdcm, [125](#)
- GDCM\_JOIN
  - gdcmStaticAssert.h, [1134](#)
- GDCM\_LEGACY
  - gdcmLegacyMacro.h, [1049](#)
- GDCM\_LEGACY\_BODY
  - gdcmLegacyMacro.h, [1049](#)
- GDCM\_LEGACY\_REPLACED\_BODY
  - gdcmLegacyMacro.h, [1049](#)
- GDCM\_LESS
  - gdcm, [125](#)
- GDCM\_LESSOREQUAL
  - gdcm, [125](#)
- GDCM\_STATIC\_ASSERT
  - gdcm::Attribute, [171](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [177](#), [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [184](#)
  - gdcmStaticAssert.h, [1134](#)
- GDCMMACROENTRY\_H
  - gdcmMacroEntry.h, [1054](#)
- GEMS
  - gdcm::Dicts, [323](#)
- GEPrivate3DModelStorage
  - gdcm::MediaStorage, [509](#)
- GRAY
  - gdcm::LookupTable, [498](#)
- GREEN
  - gdcm::LookupTable, [498](#)
- gdcm, [109](#)
  - AComp, [124](#)
  - ASComp, [124](#)
  - BOOL\_FUNCTION\_PFILE\_PFILE\_POINTER, [124](#)
  - backslash, [126](#)
  - CSComp, [124](#)
  - CompOperators, [124](#)
  - DAComp, [124](#)
  - DTCComp, [124](#)
  - eArabic, [125](#)
  - ECharSet, [125](#)
  - eCyrillic, [125](#)
  - eFind, [125](#)
  - eGB18030, [125](#)
  - eGreek, [125](#)
  - eHebrew, [125](#)
  - eImage, [125](#)
  - eJapanese, [125](#)
  - eJapaneseKanjiMultibyte, [125](#)
  - eJapaneseSupplementaryKanjiMultibyte, [125](#)
  - eKoreanHangulHanjaMultibyte, [125](#)
  - eLatin1, [125](#)
  - eLatin2, [125](#)
  - eLatin3, [125](#)
  - eLatin4, [125](#)



- eLatin5, [125](#)
- eMove, [125](#)
- ePatient, [125](#)
- ePatientRootType, [126](#)
- EQueryLevel, [125](#)
- EQueryType, [125](#)
- ERootType, [125](#)
- eSeries, [125](#)
- eStudy, [125](#)
- eStudyRootType, [126](#)
- eThai, [125](#)
- eUTF8, [125](#)
- FileList, [124](#)
- GDCM\_DIFFERENT, [124](#)
- GDCM\_EQUAL, [124](#)
- GDCM\_GREATER, [125](#)
- GDCM\_GREATEROREQUAL, [125](#)
- GDCM\_LESS, [125](#)
- GDCM\_LESSEOREQUAL, [125](#)
- GetVRFromTag, [126](#)
- GlobalInstance, [131](#)
- IconImage, [124](#)
- LD\_ALL, [126](#)
- LD\_NOSEQ, [126](#)
- LD\_NOSHADOW, [126](#)
- LD\_NOSHADOWSEQ, [126](#)
- LOComp, [124](#)
- LTComp, [124](#)
- LodModeType, [126](#)
- MacroEntry, [124](#)
- NestedMacroEntries, [124](#)
- operator!=, [126](#)
- operator<<, [126–130](#)
- operator>>, [130](#)
- operator==, [130](#)
- PNComp, [124](#)
- SHComp, [124](#)
- STComp, [124](#)
- TMComp, [124](#)
- TYPETOENCODING, [130](#)
- to\_string, [130](#)
- UIComp, [124](#)
- UTComp, [124](#)
- VRBINARY, [131](#)
- gdcmm2pnm.dox, [929](#)
- gdcmm2vtk.dox, [929](#)
- gdcmm::ASN1, [167](#)
  - ~ASN1, [168](#)
  - ASN1, [168](#)
  - ParseDump, [168](#)
  - ParseDumpFile, [168](#)
  - TestPBKDF2, [168](#)
- gdcmm::AbortEvent, [149](#)
- gdcmm::AnonymizeEvent, [151](#)
  - ~AnonymizeEvent, [153](#)
  - AnonymizeEvent, [153](#)
  - CheckEvent, [153](#)
  - GetEventName, [153](#)
  - GetTag, [153](#)
  - MakeObject, [153](#)
  - Self, [153](#)
  - SetTag, [153](#)
  - Superclass, [153](#)
- gdcmm::Anonymizer, [154](#)
  - ~Anonymizer, [156](#)
  - Anonymizer, [156](#)
  - BALCPPProtect, [156](#)
  - BasicApplicationLevelConfidentialityProfile, [156](#)
  - CanEmptyTag, [156](#)
  - ClearInternalUIDs, [157](#)
  - Empty, [157](#)
  - GetBasicApplicationLevelConfidentialityProfile←  
Attributes, [157](#)
  - GetCryptographicMessageSyntax, [157](#)
  - GetFile, [157](#)
  - New, [157](#)
  - RecurseDataSet, [157](#)
  - Remove, [157](#)
  - RemoveGroupLength, [157](#)
  - RemovePrivateTags, [157](#)
  - RemoveRetired, [158](#)
  - Replace, [158](#)
  - SetCryptographicMessageSyntax, [158](#)
  - SetFile, [158](#)
- gdcmm::AnyEvent, [158](#)
- gdcmm::ApplicationEntity, [161](#)
  - Internal, [162](#)
  - IsValid, [162](#)
  - MaxLength, [162](#)
  - MaxNumberOfComponents, [162](#)
  - Padding, [162](#)
  - Print, [162](#)
  - Separator, [162](#)
  - SetBlob, [162](#)
  - Squeeze, [162](#)
- gdcmm::Attribute
  - ArrayType, [171](#)
  - GDCM\_STATIC\_ASSERT, [171](#)
  - GetAsDataElement, [171](#)
  - GetDictVM, [171](#)
  - GetDictVR, [171](#)
  - GetNumberOfValues, [172](#)
  - GetTag, [172](#)
  - GetVM, [173](#)
  - GetVR, [173](#)
  - GetValue, [172](#)
  - GetValues, [172](#)
  - Internal, [175](#)



- operator!=, [173](#)
- operator<, [173](#)
- operator==, [173](#)
- operator[], [173](#), [174](#)
- Print, [174](#)
- Set, [174](#)
- SetByteValue, [174](#)
- SetByteValueNoSwap, [174](#)
- SetFromDataElement, [174](#)
- SetFromDataSet, [175](#)
- SetValue, [175](#)
- SetValues, [175](#)
- VMType, [171](#)
- gdcmm::Attribute< Group, Element, TVR, TVM >, [169](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, [176](#)
  - ArrayType, [177](#)
  - GDCM\_STATIC\_ASSERT, [177](#), [178](#)
  - GetAsDataElement, [178](#)
  - GetDictVM, [178](#)
  - GetDictVR, [178](#)
  - GetNumberOfValues, [178](#)
  - GetTag, [178](#)
  - GetVM, [178](#)
  - GetVR, [178](#)
  - GetValue, [178](#)
  - GetValues, [178](#)
  - Internal, [180](#)
  - operator!=, [179](#)
  - operator<, [179](#)
  - operator==, [179](#)
  - Print, [179](#)
  - Set, [179](#)
  - SetByteValue, [179](#)
  - SetByteValueNoSwap, [179](#)
  - SetFromDataElement, [179](#)
  - SetFromDataSet, [180](#)
  - SetValue, [180](#)
  - VMType, [177](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM1\_3 >, [180](#)
  - GetVM, [181](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM1\_8 >, [181](#)
  - GetVM, [182](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, [183](#)
  - ~Attribute, [184](#)
  - ArrayType, [184](#)
  - Attribute, [184](#)
  - GDCM\_STATIC\_ASSERT, [184](#)
  - GetAsDataElement, [184](#)
  - GetDictVM, [184](#)
  - GetDictVR, [185](#)
  - GetNumberOfValues, [185](#)
  - GetTag, [185](#)
  - GetVM, [185](#)
  - GetVR, [185](#)
  - GetValue, [185](#)
  - GetValues, [185](#)
  - operator[], [185](#)
  - Print, [185](#)
  - Set, [186](#)
  - SetByteValue, [186](#)
  - SetFromDataElement, [186](#)
  - SetFromDataSet, [186](#)
  - SetNumberOfValues, [186](#)
  - SetValue, [186](#)
  - SetValues, [187](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM2\_2n >, [187](#)
  - GetVM, [188](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM2\_n >, [188](#)
  - GetVM, [189](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM3\_3n >, [190](#)
  - GetVM, [191](#)
- gdcmm::Attribute< Group, Element, TVR, VM::VM3\_n >, [191](#)
  - GetVM, [192](#)
- gdcmm::AudioCodec, [193](#)
  - ~AudioCodec, [194](#)
  - AudioCodec, [194](#)
  - CanCode, [194](#)
  - CanDecode, [195](#)
  - Decode, [195](#)
- gdcmm::Base64, [195](#)
  - Decode, [195](#)
  - Encode, [196](#)
  - GetDecodeLength, [196](#)
  - GetEncodeLength, [196](#)
- gdcmm::BaseRootQuery, [200](#)
  - ~BaseRootQuery, [202](#)
  - AddQueryDataSet, [202](#)
  - BaseRootQuery, [202](#)
  - Construct, [203](#)
  - GetAbstractSyntaxUID, [203](#)
  - GetQueryDataSet, [203](#)
  - GetQueryLevelFromQueryRoot, [203](#)
  - GetQueryLevelFromString, [203](#)
  - GetQueryLevelString, [203](#)
  - GetTagListByLevel, [203](#)
  - InitializeDataSet, [203](#)
  - mDataSet, [204](#)
  - mHelpDescription, [204](#)
  - mImage, [204](#)
  - mPatient, [204](#)
  - mRootType, [204](#)

- mSeries, [204](#)
- mStudy, [204](#)
- Print, [203](#)
- QueryFactory, [204](#)
- SetSearchParameter, [203](#)
- ValidateQuery, [203](#)
- WriteHelpFile, [204](#)
- WriteQuery, [204](#)
- gdcm::BasicOffsetTable, [207](#)
  - BasicOffsetTable, [208](#)
  - operator<<, [209](#)
  - Read, [208](#)
- gdcm::Bitmap, [209](#)
  - ~Bitmap, [212](#)
  - AreOverlaysInPixelData, [212](#)
  - Bitmap, [212](#)
  - Clear, [212](#)
  - ComputeLossyFlag, [212](#)
  - Dimensions, [216](#)
  - GetBuffer, [212](#)
  - GetBuffer2, [212](#)
  - GetBufferLength, [212](#)
  - GetColumns, [213](#)
  - GetDataElement, [213](#)
  - GetDimension, [213](#)
  - GetDimensions, [213](#)
  - GetLUT, [213](#)
  - GetNeedByteSwap, [213](#)
  - GetNumberOfDimensions, [213](#)
  - GetPhotometricInterpretation, [213](#)
  - GetPixelFormat, [213, 214](#)
  - GetPlanarConfiguration, [214](#)
  - GetRows, [214](#)
  - GetTransferSyntax, [214](#)
  - ImageChangeTransferSyntax, [216](#)
  - IsEmpty, [214](#)
  - IsLossy, [214](#)
  - IsTransferSyntaxCompatible, [214](#)
  - LUT, [216](#)
  - LUTPtr, [212](#)
  - LossyFlag, [216](#)
  - NeedByteSwap, [216](#)
  - NumberOfDimensions, [216](#)
  - PF, [216](#)
  - PI, [216](#)
  - PixelData, [216](#)
  - PixmapReader, [216](#)
  - PlanarConfiguration, [216](#)
  - Print, [214](#)
  - SetColumns, [214](#)
  - SetDataElement, [214](#)
  - SetDimension, [214](#)
  - SetDimensions, [215](#)
  - SetLUT, [215](#)
  - SetLossyFlag, [215](#)
  - SetNeedByteSwap, [215](#)
  - SetNumberOfDimensions, [215](#)
  - SetPhotometricInterpretation, [215](#)
  - SetPixelFormat, [215](#)
  - SetPlanarConfiguration, [215](#)
  - SetRows, [215](#)
  - SetTransferSyntax, [216](#)
  - TS, [217](#)
  - TryJPEG2000Codec, [216](#)
  - TryJPEG2000Codec2, [216](#)
  - TryJPEGCodec, [216](#)
  - TryJPEGCodec2, [216](#)
  - TryJPEGLSCodec, [216](#)
  - TryKAKADUCodec, [216](#)
  - TryPVRGCodec, [216](#)
  - TryRAWCodec, [216](#)
  - TryRLECodec, [216](#)
- gdcm::BitmapToBitmapFilter, [217](#)
  - ~BitmapToBitmapFilter, [218](#)
  - BitmapToBitmapFilter, [218](#)
  - GetOutput, [218](#)
  - GetOutputAsBitmap, [218](#)
  - Input, [218](#)
  - Output, [218](#)
  - SetInput, [218](#)
- gdcm::BoxRegion, [219](#)
  - ~BoxRegion, [220](#)
  - Area, [220](#)
  - BoundingBox, [220](#)
  - BoxRegion, [220](#)
  - Clone, [221](#)
  - ComputeBoundingBox, [221](#)
  - Empty, [221](#)
  - GetXMax, [221](#)
  - GetXMin, [221](#)
  - GetYMax, [221](#)
  - GetYMin, [221](#)
  - GetZMax, [221](#)
  - GetZMin, [221](#)
  - IsValid, [221](#)
  - operator=, [221](#)
  - Print, [221](#)
  - SetDomain, [221](#)
- gdcm::ByteBuffer, [222](#)
  - ByteBuffer, [222](#)
  - Get, [222](#)
  - GetStart, [222](#)
  - ShiftEnd, [222](#)
  - UpdatePosition, [222](#)
- gdcm::ByteSwap
  - Swap, [223](#)
  - SwapFromSwapCodeIntoSystem, [223](#)
  - SwapRange, [223](#)

- SwapRangeFromSwapCodeIntoSystem, [223](#)
- SystemIsBigEndian, [223](#)
- SystemIsLittleEndian, [223](#)
- gdcmm::ByteSwap< T >, [223](#)
- gdcmm::ByteSwapFilter, [224](#)
  - ~ByteSwapFilter, [224](#)
  - ByteSwap, [224](#)
  - ByteSwapFilter, [224](#)
  - SetByteSwapTag, [224](#)
- gdcmm::ByteValue, [224](#)
  - ~ByteValue, [227](#)
  - Append, [227](#)
  - ByteValue, [226](#)
  - Clear, [227](#)
  - ComputeLength, [227](#)
  - Fill, [227](#)
  - GetBuffer, [227](#)
  - GetLength, [227](#)
  - GetPointer, [227](#)
  - IsEmpty, [228](#)
  - IsPrintable, [228](#)
  - operator const std::vector< char > &, [228](#)
  - operator=, [228](#)
  - operator==, [228](#)
  - Print, [228](#)
  - PrintASCII, [228](#)
  - PrintASCIIXML, [228](#)
  - PrintGroupLength, [228](#)
  - PrintHex, [228](#)
  - PrintHexXML, [228](#)
  - PrintPNXML, [228](#)
  - Read, [228](#), [229](#)
  - SetLength, [229](#)
  - SetLengthOnly, [229](#)
  - Write, [229](#)
  - WriteBuffer, [229](#)
- gdcmm::CAPICryptoFactory, [229](#)
  - CAPICryptoFactory, [230](#)
  - CreateCMSProvider, [230](#)
- gdcmm::CAPICryptographicMessageSyntax, [230](#)
  - ~CAPICryptographicMessageSyntax, [232](#)
  - CAPICryptographicMessageSyntax, [232](#)
  - Decrypt, [232](#)
  - Encrypt, [232](#)
  - GetCipherType, [232](#)
  - GetInitialized, [232](#)
  - ParseCertificateFile, [232](#)
  - ParseKeyFile, [232](#)
  - SetCipherType, [232](#)
  - SetPassword, [232](#)
- gdcmm::CP246ExplicitDataElement, [258](#)
  - GetLength, [259](#)
  - Read, [259](#)
  - ReadPreValue, [259](#)
  - ReadValue, [260](#)
  - ReadWithLength, [260](#)
- gdcmm::CSAElement, [264](#)
  - CSAElement, [266](#)
  - DataField, [268](#)
  - DataPtr, [266](#)
  - GetByteValue, [266](#)
  - GetKey, [266](#)
  - GetName, [266](#)
  - GetNoOfItems, [266](#)
  - GetSyngoDT, [266](#)
  - GetVM, [267](#)
  - GetVR, [267](#)
  - GetValue, [266](#)
  - IsEmpty, [267](#)
  - KeyField, [268](#)
  - NameField, [268](#)
  - NoOfItemsField, [268](#)
  - operator<, [267](#)
  - operator<<, [268](#)
  - operator=, [267](#)
  - operator==, [267](#)
  - SetByteValue, [267](#)
  - SetKey, [267](#)
  - SetName, [267](#)
  - SetNoOfItems, [267](#)
  - SetSyngoDT, [267](#)
  - SetVM, [267](#)
  - SetVR, [267](#)
  - SetValue, [267](#)
  - SyngoDTField, [268](#)
  - VRField, [268](#)
  - ValueMultiplicityField, [268](#)
- gdcmm::CSAHeader, [268](#)
  - ~CSAHeader, [270](#)
  - CSAHeader, [270](#)
  - CSAHeaderType, [270](#)
  - DATASET\_FORMAT, [270](#)
  - FindCSAElementByName, [270](#)
  - GetCSADataInfo, [271](#)
  - GetCSAEEnd, [271](#)
  - GetCSAElementByName, [271](#)
  - GetCSAImageHeaderInfoTag, [271](#)
  - GetCSASeriesHeaderInfoTag, [271](#)
  - GetDataSet, [271](#)
  - GetFormat, [272](#)
  - GetInterfile, [272](#)
  - INTERFILE, [270](#)
  - LoadFromDataElement, [272](#)
  - NOMAGIC, [270](#)
  - operator<<, [272](#)
  - Print, [272](#)
  - Read, [272](#)
  - SV10, [270](#)

- UNKNOWN, [270](#)
- Write, [272](#)
- ZEROED\_OUT, [270](#)
- gdcmm::CSAHeaderDict, [272](#)
  - AddCSAHeaderDictEntry, [273](#)
  - Begin, [273](#)
  - CSAHeaderDict, [273](#)
  - ConstIterator, [273](#)
  - Dicts, [274](#)
  - End, [274](#)
  - GetCSAHeaderDictEntry, [274](#)
  - IsEmpty, [274](#)
  - Iterator, [273](#)
  - LoadDefault, [274](#)
  - MapCSAHeaderDictEntry, [273](#)
  - operator<<, [274](#)
- gdcmm::CSAHeaderDictEntry, [274](#)
  - CSAHeaderDictEntry, [275](#)
  - GetDescription, [275](#)
  - GetName, [275](#)
  - GetVM, [275](#)
  - GetVR, [275](#)
  - operator<, [275](#)
  - operator<<, [276](#)
  - SetDescription, [275](#)
  - SetName, [276](#)
  - SetVM, [276](#)
  - SetVR, [276](#)
- gdcmm::CSAHeaderDictException, [276](#)
- gdcmm::CodeString, [245](#)
  - CodeString, [247](#)
  - const\_iterator, [247](#)
  - const\_reference, [247](#)
  - const\_reverse\_iterator, [247](#)
  - difference\_type, [247](#)
  - GetAsString, [247](#)
  - IsValid, [248](#)
  - iterator, [247](#)
  - operator!=, [248](#)
  - operator<<, [248](#)
  - operator==, [248](#)
  - pointer, [247](#)
  - reference, [247](#)
  - reverse\_iterator, [247](#)
  - Size, [248](#)
  - size\_type, [247](#)
  - TrimInternal, [248](#)
  - value\_type, [247](#)
- gdcmm::Codec, [243](#)
- gdcmm::Coder, [244](#)
  - ~Coder, [245](#)
  - CanCode, [245](#)
  - Code, [245](#)
  - InternalCode, [245](#)
- gdcmm::Command, [248](#)
  - ~Command, [250](#)
  - Command, [250](#)
  - Execute, [250](#)
- gdcmm::CommandDataSet, [250](#)
  - ~CommandDataSet, [252](#)
  - CommandDataSet, [252](#)
  - Insert, [252](#)
  - operator<<, [252](#)
  - Read, [252](#)
  - Replace, [252](#)
  - Write, [252](#)
- gdcmm::CompositeNetworkFunctions, [253](#)
  - CEcho, [254](#)
  - CFind, [256](#)
  - CMove, [256](#)
  - CStore, [257](#)
  - ConstructQuery, [257](#)
  - KeyValuePairArrayType, [254](#)
  - KeyValuePairType, [254](#)
- gdcmm::ConstCharWrapper, [258](#)
  - ConstCharWrapper, [258](#)
  - operator const char \*, [258](#)
- gdcmm::CryptoFactory, [260](#)
  - ~CryptoFactory, [261](#)
  - CAPI, [261](#)
  - CreateCMSProvider, [261](#)
  - CryptoFactory, [261](#)
  - CryptoLib, [261](#)
  - DEFAULT, [261](#)
  - GetFactoryInstance, [261](#)
  - OPENSSL, [261](#)
  - OPENSSL7, [261](#)
- gdcmm::CryptographicMessageSyntax, [262](#)
  - ~CryptographicMessageSyntax, [263](#)
  - AES128\_CIPHER, [262](#)
  - AES192\_CIPHER, [262](#)
  - AES256\_CIPHER, [263](#)
  - CipherTypes, [262](#)
  - CryptographicMessageSyntax, [263](#)
  - DES3\_CIPHER, [262](#)
  - Decrypt, [263](#)
  - Encrypt, [263](#)
  - GetCipherType, [263](#)
  - ParseCertificateFile, [263](#)
  - ParseKeyFile, [263](#)
  - SetCipherType, [263](#)
  - SetPassword, [263](#)
- gdcmm::Curve, [280](#)
  - ~Curve, [281](#)
  - Curve, [281](#)
  - Decode, [281](#)
  - GetAsPoints, [281](#)
  - GetCurveDataDescriptor, [281](#)

- GetDataValueRepresentation, [282](#)
- GetDimensions, [282](#)
- GetGroup, [282](#)
- GetNumberOfCurves, [282](#)
- GetNumberOfPoints, [282](#)
- GetTypeOfData, [282](#)
- GetTypeOfDataDescription, [282](#)
- IsEmpty, [282](#)
- Print, [282](#)
- SetCoordinateStartValue, [282](#)
- SetCoordinateStepValue, [282](#)
- SetCurve, [282](#)
- SetCurveDataDescriptor, [282](#)
- SetCurveDescription, [282](#)
- SetDataValueRepresentation, [282](#)
- SetDimensions, [282](#)
- SetGroup, [282](#)
- SetNumberOfPoints, [282](#)
- SetTypeOfData, [282](#)
- Update, [282](#)
- gdcm::DICOMDIR, [310](#)
  - DICOMDIR, [311](#)
- gdcm::DICOMDIRGenerator, [311](#)
  - ~DICOMDIRGenerator, [312](#)
  - AddImageDirectoryRecord, [312](#)
  - AddPatientDirectoryRecord, [312](#)
  - AddSeriesDirectoryRecord, [312](#)
  - AddStudyDirectoryRecord, [312](#)
  - DICOMDIRGenerator, [312](#)
  - FilenameType, [312](#)
  - FileNamesType, [312](#)
  - Generate, [312](#)
  - GetFile, [312](#)
  - GetScanner, [313](#)
  - SetDescriptor, [313](#)
  - SetFile, [313](#)
  - SetFileNames, [313](#)
  - SetRootDirectory, [313](#)
- gdcm::DataElement, [282](#)
  - Clear, [286](#)
  - DataElement, [286](#)
  - Empty, [286](#)
  - GetByteValue, [286](#)
  - GetLength, [286](#)
  - GetSequenceOfFragments, [286](#)
  - GetTag, [287](#)
  - GetVL, [287](#)
  - GetVR, [288](#)
  - GetValue, [287](#)
  - GetValueAsSQ, [287](#)
  - IsEmpty, [288](#)
  - IsUndefinedLength, [288](#)
  - operator<, [288](#)
  - operator<<, [291](#)
  - operator=, [288](#)
  - operator==, [288](#)
  - Read, [288](#)
  - ReadOrSkip, [289](#)
  - ReadPreValue, [289](#)
  - ReadValue, [289](#)
  - ReadValueWithLength, [289](#)
  - ReadWithLength, [289](#)
  - SetByteValue, [289](#)
  - SetTag, [289](#)
  - SetVL, [290](#)
  - SetVLToUndefined, [290](#)
  - SetVR, [290](#)
  - SetValue, [289](#)
  - SetValueFieldLength, [290](#)
  - TagField, [291](#)
  - VRField, [291](#)
  - ValueField, [291](#)
  - ValueLengthField, [291](#)
  - ValuePtr, [285](#)
  - Write, [290](#)
- gdcm::DataElementException, [291](#)
- gdcm::DataEvent, [292](#)
  - ~DataEvent, [294](#)
  - CheckEvent, [294](#)
  - DataEvent, [293](#), [294](#)
  - GetData, [294](#)
  - GetDataLength, [294](#)
  - GetEventName, [294](#)
  - MakeObject, [294](#)
  - Self, [293](#)
  - SetData, [294](#)
  - Superclass, [293](#)
- gdcm::DataSet, [294](#)
  - Begin, [297](#)
  - CSAHeader, [301](#)
  - Clear, [297](#)
  - ComputeDataElement, [297](#)
  - ComputeGroupLength, [298](#)
  - ConstIterator, [297](#)
  - DataElementSet, [297](#)
  - End, [298](#)
  - FindDataElement, [298](#)
  - FindNextDataElement, [298](#)
  - GetDEEnd, [299](#)
  - GetDES, [299](#)
  - GetDataElement, [298](#), [299](#)
  - GetLength, [299](#)
  - GetMediaStorage, [299](#)
  - GetPrivateCreator, [299](#)
  - Insert, [299](#)
  - InsertDataElement, [299](#)
  - IsEmpty, [299](#)
  - Iterator, [297](#)

- operator<<, 301
- operator(), 300
- operator=, 300
- operator[], 300
- Print, 300
- Read, 300
- ReadNested, 300
- ReadSelectedPrivateTags, 300
- ReadSelectedPrivateTagsWithLength, 300
- ReadSelectedTags, 300
- ReadSelectedTagsWithLength, 300
- ReadUpToTag, 300
- ReadUpToTagWithLength, 300
- ReadWithLength, 300
- Remove, 300
- Replace, 300
- ReplaceEmpty, 301
- Size, 301
- SizeType, 297
- Write, 301
- gdcmm::DataSetEvent, 301
  - ~DataSetEvent, 303
  - CheckEvent, 303
  - DataSetEvent, 303
  - GetDataSet, 303
  - GetEventName, 303
  - MakeObject, 303
  - Self, 303
  - Superclass, 303
- gdcmm::DataSetHelper, 304
  - ComputeVR, 304
- gdcmm::Decoder, 304
  - ~Decoder, 305
  - CanDecode, 305
  - Decode, 305
  - DecodeByStreams, 305
- gdcmm::DefinedTerms, 306
  - DefinedTerms, 306
- gdcmm::Defs, 306
  - ~Defs, 307
  - Defs, 307
  - GetIODFromFile, 307
  - GetIODNameFromMediaStorage, 307
  - GetIODs, 307, 308
  - GetMacros, 308
  - GetModules, 308
  - GetTypeFromTag, 308
  - Global, 308
  - IsEmpty, 308
  - LoadDefaults, 308
  - LoadFromFile, 308
  - Verify, 308
- gdcmm::DeltaEncodingCodec, 309
  - ~DeltaEncodingCodec, 310
- CanDecode, 310
- Decode, 310
- DeltaEncodingCodec, 310
- gdcmm::Dict, 313
  - AddDictEntry, 314
  - Begin, 314
  - ConstIterator, 314
  - Dict, 314
  - Dicts, 315
  - End, 314
  - GetDictEntry, 315
  - GetDictEntryByKeyword, 315
  - GetDictEntryByName, 315
  - GetKeywordFromTag, 315
  - IsEmpty, 315
  - Iterator, 314
  - LoadDefault, 315
  - MapDictEntry, 314
  - operator<<, 315
- gdcmm::DictConverter, 316
  - ~DictConverter, 317
  - AddGroupLength, 317
  - Convert, 317
  - ConvertToCXX, 317
  - ConvertToXML, 317
  - DICT\_DEBUG, 317
  - DICT\_DEFAULT, 317
  - DICT\_XML, 317
  - DictConverter, 317
  - GetDictName, 317
  - GetInputFilename, 317
  - GetOutputFilename, 317
  - GetOutputType, 317
  - OutputTypes, 317
  - ReadVM, 317
  - ReadVR, 317
  - Readuint16, 317
  - SetDictName, 317
  - SetInputFileName, 317
  - SetOutputFileName, 317
  - SetOutputType, 317
  - WriteFooter, 317
  - WriteHeader, 318
- gdcmm::DictEntry, 318
  - Dict, 320
  - DictEntry, 319
  - GetKeyword, 319
  - GetName, 319
  - GetRetired, 319
  - GetVM, 319
  - GetVR, 319
  - IsUnique, 320
  - operator<<, 320
  - SetElementXX, 320

- SetGroupXX, [320](#)
- SetKeyword, [320](#)
- SetName, [320](#)
- SetRetired, [320](#)
- SetVM, [320](#)
- SetVR, [320](#)
- gdcm::DictPrinter, [320](#)
  - ~DictPrinter, [322](#)
  - DictPrinter, [322](#)
  - Print, [322](#)
  - PrintDataElement2, [322](#)
  - PrintDataSet2, [322](#)
- gdcm::Dicts, [322](#)
  - ~Dicts, [323](#)
  - ConstructorType, [323](#)
  - Dicts, [323](#)
  - GEMS, [323](#)
  - GetCSAHeaderDict, [324](#)
  - GetConstructorString, [323](#)
  - GetDictEntry, [324](#)
  - GetPrivateDict, [324](#)
  - GetPublicDict, [324](#)
  - Global, [324](#)
  - IsEmpty, [324](#)
  - LoadDefaults, [324](#)
  - operator<<, [324](#)
  - PHILIPS, [323](#)
  - SIEMENS, [323](#)
- gdcm::DirectionCosines, [326](#)
  - ~DirectionCosines, [327](#)
  - ComputeDistAlongNormal, [327](#)
  - Cross, [327](#)
  - CrossDot, [327](#)
  - DirectionCosines, [327](#)
  - Dot, [327](#)
  - IsValid, [327](#)
  - Normalize, [327](#)
  - operator const double \*, [328](#)
  - Print, [328](#)
  - SetFromString, [328](#)
- gdcm::Directory, [328](#)
  - ~Directory, [329](#)
  - Directory, [329](#)
  - Explore, [329](#)
  - FilenameType, [329](#)
  - FileNamesType, [329](#)
  - GetDirectories, [329](#)
  - GetFileNames, [330](#)
  - GetToplevel, [330](#)
  - Load, [330](#)
  - operator<<, [330](#)
  - Print, [330](#)
- gdcm::DirectoryHelper, [331](#)
  - GetCTImageSeriesUIDs, [331](#)
  - GetFileNamesFromSeriesUIDs, [331](#)
  - GetFrameOfReference, [331](#)
  - GetMRIImageSeriesUIDs, [331](#)
  - GetRTStructSeriesUIDs, [331](#)
  - GetSOPClassUID, [332](#)
  - GetSeriesUIDsBySOPClassUID, [332](#)
  - GetStringValueFromTag, [332](#)
  - LoadImageFromFiles, [332](#)
  - RetrieveSOPInstanceUIDFromIndex, [332](#)
  - RetrieveSOPInstanceUIDFromZPosition, [332](#)
- gdcm::DummyValueGenerator, [332](#)
  - Generate, [332](#)
- gdcm::Dumper, [333](#)
  - ~Dumper, [334](#)
  - Dumper, [334](#)
- gdcm::Element
  - GetAsDataElement, [337](#)
  - GetLength, [337](#)
  - GetVM, [337](#)
  - GetVR, [337](#)
  - GetValue, [337](#)
  - GetValues, [337](#)
  - Internal, [337](#)
  - operator[], [337](#)
  - Print, [337](#)
  - Read, [337](#)
  - Set, [337](#)
  - SetFromDataElement, [337](#)
  - SetNoSwap, [337](#)
  - SetValue, [337](#)
  - Type, [337](#)
  - Write, [337](#)
- gdcm::Element< TVR, TVM >, [335](#)
- gdcm::Element< TVR, VM::VM1\_2 >, [338](#)
  - Parent, [339](#)
  - SetLength, [339](#)
- gdcm::Element< TVR, VM::VM1\_n >, [339](#)
  - ~Element, [340](#)
  - Element, [340](#)
  - GetAsDataElement, [340](#)
  - GetLength, [340](#)
  - GetVM, [340](#)
  - GetVR, [340](#)
  - GetValue, [340](#)
  - operator=, [341](#)
  - operator[], [341](#)
  - Print, [341](#)
  - Read, [341](#)
  - Set, [341](#)
  - SetArray, [341](#)
  - SetFromDataElement, [341](#)
  - SetLength, [341](#)
  - SetNoSwap, [341](#)
  - SetValue, [341](#)

- Type, [340](#)
- Write, [341](#)
- WriteASCII, [341](#)
- gdcmm::Element< TVR, VM::VM2\_2n >, [341](#)
  - Parent, [343](#)
  - SetLength, [343](#)
- gdcmm::Element< TVR, VM::VM2\_n >, [343](#)
  - Parent, [344](#)
  - SetLength, [344](#)
- gdcmm::Element< TVR, VM::VM3\_3n >, [344](#)
  - Parent, [346](#)
  - SetLength, [346](#)
- gdcmm::Element< TVR, VM::VM3\_n >, [346](#)
  - Parent, [347](#)
  - SetLength, [347](#)
- gdcmm::Element< VR::AS, VM::VM5 >, [347](#)
  - GetLength, [348](#)
  - Internal, [348](#)
  - Print, [348](#)
- gdcmm::Element< VR::OB, VM::VM1 >, [348](#)
- gdcmm::Element< VR::OW, VM::VM1 >, [349](#)
- gdcmm::ElementDisableCombinations< TVR, TVM >, [351](#)
- gdcmm::ElementDisableCombinations< VR::OB, VM::V←  
M1\_n >, [352](#)
- gdcmm::ElementDisableCombinations< VR::OW, VM::V←  
M1\_n >, [352](#)
- gdcmm::EncapsulatedDocument, [352](#)
  - EncapsulatedDocument, [352](#)
- gdcmm::EncodingImplementation< T >, [353](#)
- gdcmm::EncodingImplementation< VR::VRASCII >, [353](#)
  - Read, [353](#)
  - ReadComputeLength, [353](#)
  - ReadNoSwap, [354](#)
  - Write, [354](#)
- gdcmm::EncodingImplementation< VR::VRBINARY >, [354](#)
  - Read, [354](#)
  - ReadComputeLength, [354](#)
  - ReadNoSwap, [355](#)
  - Write, [355](#)
- gdcmm::EndEvent, [355](#)
- gdcmm::EnumeratedValues, [356](#)
  - EnumeratedValues, [357](#)
- gdcmm::Event, [357](#)
  - ~Event, [358](#)
  - CheckEvent, [358](#)
  - Event, [358](#)
  - GetEventName, [358](#)
  - MakeObject, [358](#)
  - Print, [358](#)
- gdcmm::Exception, [359](#)
  - ~Exception, [360](#)
  - Exception, [360](#)
  - GetDescription, [360](#)
  - what, [360](#)
- gdcmm::ExitEvent, [361](#)
- gdcmm::ExplicitDataElement, [362](#)
  - GetLength, [363](#)
  - Read, [363](#)
  - ReadPreValue, [363](#)
  - ReadValue, [363](#)
  - ReadWithLength, [363](#)
  - Write, [363](#)
- gdcmm::ExplicitImplicitDataElement, [363](#)
  - GetLength, [365](#)
  - Read, [365](#)
  - ReadPreValue, [365](#)
  - ReadValue, [365](#)
  - ReadWithLength, [365](#)
- gdcmm::Fiducials, [365](#)
  - Fiducials, [365](#)
- gdcmm::File, [366](#)
  - File, [367](#)
  - GetDataSet, [367](#), [368](#)
  - GetHeader, [368](#)
  - operator<<, [369](#)
  - Read, [368](#)
  - SetDataSet, [368](#)
  - SetHeader, [368](#)
  - Write, [368](#)
- gdcmm::FileAnonymizer, [369](#)
  - ~FileAnonymizer, [371](#)
  - Empty, [371](#)
  - FileAnonymizer, [371](#)
  - Remove, [371](#)
  - Replace, [371](#)
  - SetInputFileName, [371](#)
  - SetOutputFileName, [371](#)
  - Write, [372](#)
- gdcmm::FileChangeTransferSyntax, [372](#)
  - ~FileChangeTransferSyntax, [374](#)
  - Change, [374](#)
  - FileChangeTransferSyntax, [374](#)
  - GetCodec, [374](#)
  - New, [374](#)
  - SetInputFileName, [374](#)
  - SetOutputFileName, [374](#)
  - SetTransferSyntax, [374](#)
- gdcmm::FileDerivation, [374](#)
  - ~FileDerivation, [375](#)
  - AddDerivationDescription, [375](#)
  - AddPurposeOfReferenceCodeSequence, [375](#)
  - AddReference, [375](#)
  - AddSourceImageSequence, [376](#)
  - Derive, [376](#)
  - FileDerivation, [375](#)
  - GetFile, [376](#)
  - SetDerivationCodeSequenceCodeValue, [376](#)
  - SetDerivationDescription, [376](#)



- SetFile, [376](#)
- SetPurposeOfReferenceCodeSequenceCodeValue, [376](#)
- gdcm::FileExplicitFilter, [377](#)
  - ~FileExplicitFilter, [378](#)
  - Change, [378](#)
  - ChangeFMI, [378](#)
  - FileExplicitFilter, [378](#)
  - GetFile, [378](#)
  - ProcessDataSet, [378](#)
  - SetChangePrivateTags, [378](#)
  - SetFile, [378](#)
  - SetRecomputeItemLength, [378](#)
  - SetRecomputeSequenceLength, [378](#)
  - SetUseVRUN, [379](#)
- gdcm::FileMetaInformation, [379](#)
  - AppendImplementationClassUID, [381](#)
  - ComputeDataSetMediaStorageSOPClass, [381](#)
  - ComputeDataSetTransferSyntax, [381](#)
  - DataSetMS, [383](#)
  - DataSetTS, [383](#)
  - Default, [381](#)
  - FileMetaInformation, [381](#)
  - FillFromDataSet, [381](#)
  - GetDataSetTransferSyntax, [381](#)
  - GetFileMetaInformationVersion, [382](#)
  - GetFullLength, [382](#)
  - GetGDCMImplementationClassUID, [382](#)
  - GetGDCMImplementationVersionName, [382](#)
  - GetGDCMSourceApplicationEntityTitle, [382](#)
  - GetImplementationClassUID, [382](#)
  - GetImplementationVersionName, [382](#)
  - GetMediaStorage, [382](#)
  - GetMediaStorageAsString, [382](#)
  - GetMetaInformationTS, [382](#)
  - GetPreamble, [382](#)
  - GetSourceApplicationEntityTitle, [382](#)
  - Insert, [382](#)
  - IsValid, [382](#)
  - MetaInformationTS, [384](#)
  - operator<<, [383](#)
  - Read, [382](#)
  - ReadCompat, [382](#)
  - ReadCompatInternal, [383](#)
  - Replace, [383](#)
  - SetDataSetTransferSyntax, [383](#)
  - SetImplementationClassUID, [383](#)
  - SetImplementationVersionName, [383](#)
  - SetPreamble, [383](#)
  - SetSourceApplicationEntityTitle, [383](#)
  - Write, [383](#)
- gdcm::FileNameEvent, [386](#)
  - ~FileNameEvent, [388](#)
  - CheckEvent, [388](#)
  - FileNameEvent, [388](#)
  - GetEventName, [388](#)
  - GetFileName, [388](#)
  - MakeObject, [388](#)
  - Self, [388](#)
  - SetFileName, [388](#)
  - Superclass, [388](#)
- gdcm::FileSet, [391](#)
  - AddFile, [391](#), [392](#)
  - FileSet, [391](#)
  - FileType, [391](#)
  - FilesType, [391](#)
  - GetFiles, [392](#)
  - operator<<, [392](#)
  - SetFiles, [392](#)
- gdcm::FileStreamer, [392](#)
  - ~FileStreamer, [394](#)
  - AppendToDataElement, [394](#)
  - AppendToGroupDataElement, [394](#)
  - CheckDataElement, [394](#)
  - CheckTemplateFileName, [394](#)
  - FileStreamer, [394](#)
  - New, [394](#)
  - ReserveDataElement, [394](#)
  - ReserveGroupDataElement, [395](#)
  - SetOutputFileName, [395](#)
  - SetTemplateFileName, [395](#)
  - StartDataElement, [395](#)
  - StartGroupDataElement, [395](#)
  - StopDataElement, [395](#)
  - StopGroupDataElement, [395](#)
- gdcm::FileWithName, [396](#)
  - FileWithName, [397](#)
  - filename, [397](#)
- gdcm::Filename, [384](#)
  - EndWith, [385](#)
  - Filename, [385](#)
  - GetExtension, [385](#)
  - GetFileName, [385](#)
  - GetName, [385](#)
  - GetPath, [385](#)
  - IsEmpty, [385](#)
  - IsIdentical, [385](#)
  - Join, [385](#)
  - operator const char \*, [386](#)
  - ToUnixSlashes, [386](#)
  - ToWindowsSlashes, [386](#)
- gdcm::FilenameGenerator, [388](#)
  - ~FilenameGenerator, [390](#)
  - FilenameGenerator, [389](#)
  - FilenameType, [389](#)
  - FilenamesType, [389](#)
  - Generate, [390](#)
  - GetFilename, [390](#)

- GetFileNames, 390
- GetNumberOfFileNames, 390
- GetPattern, 390
- GetPrefix, 390
- SetNumberOfFileNames, 390
- SetPattern, 390
- SetPrefix, 390
- SizeType, 389
- gdcmm::FindPatientRootQuery, 397
  - FindPatientRootQuery, 398
  - GetAbstractSyntaxUID, 398
  - GetTagListByLevel, 398
  - InitializeDataSet, 399
  - QueryFactory, 399
  - ValidateQuery, 399
- gdcmm::FindStudyRootQuery, 399
  - FindStudyRootQuery, 401
  - GetAbstractSyntaxUID, 401
  - GetTagListByLevel, 401
  - InitializeDataSet, 401
  - QueryFactory, 401
  - ValidateQuery, 401
- gdcmm::Fragment, 401
  - ComputeLength, 403
  - Fragment, 403
  - GetLength, 403
  - operator<<, 404
  - Read, 403
  - ReadBacktrack, 403
  - ReadPreValue, 403
  - ReadValue, 403
  - Write, 403
- gdcmm::Global, 404
  - ~Global, 405
  - Append, 405
  - GetDefs, 405
  - GetDicts, 405
  - GetInstance, 405
  - Global, 405
  - LoadResourcesFiles, 406
  - Locate, 406
  - operator<<, 406
  - Prepend, 406
- gdcmm::GroupDict, 406
  - ~GroupDict, 407
  - Add, 407
  - GetAbbreviation, 407
  - GetName, 407
  - GroupDict, 407
  - GroupStringVector, 407
  - Insert, 407
  - operator<<, 408
  - Size, 408
- gdcmm::IOD, 459
  - AddIODEntry, 460
  - Clear, 460
  - GetIODEntry, 460
  - GetNumberOfIODs, 460
  - GetTypeFromTag, 460
  - IOD, 459
  - MapIODEntry, 459
  - operator<<, 460
  - SizeType, 459
- gdcmm::IODEntry, 460
  - GetIE, 461
  - GetName, 461
  - GetRef, 461
  - GetUsage, 461
  - GetUsageType, 462
  - IODEntry, 461
  - operator<<, 462
  - SetIE, 462
  - SetName, 462
  - SetRef, 462
  - SetUsage, 462
- gdcmm::IODs, 462
  - AddIOD, 463
  - Begin, 463
  - Clear, 463
  - End, 463
  - GetIOD, 463
  - IODMapType, 463
  - IODMapTypeConstIterator, 463
  - IODName, 463
  - IODs, 463
  - operator<<, 463
- gdcmm::IPPSorter, 464
  - ComputeZSpacing, 467
  - DirCosTolerance, 467
  - DropDuplicatePositions, 467
  - GetDirectionCosinesTolerance, 465
  - GetZSpacing, 465
  - GetZSpacingTolerance, 466
  - IPPSorter, 465
  - SetComputeZSpacing, 466
  - SetDirectionCosinesTolerance, 466
  - SetDropDuplicatePositions, 466
  - SetZSpacingTolerance, 466
  - Sort, 466
  - ZSpacing, 467
  - ZTolerance, 467
- gdcmm::IconImageFilter, 408
  - ~IconImageFilter, 409
  - Extract, 409
  - ExtractIconImages, 409
  - ExtractVeprolIconImages, 409
  - GetFile, 409
  - GetIconImage, 409

- GetNumberOfIconImages, 410
- IconImageFilter, 409
- SetFile, 410
- gdcmm::IconImageGenerator, 410
  - ~IconImageGenerator, 411
  - AutoPixelMinMax, 411
  - ConvertRGBToPaletteColor, 411
  - Generate, 411
  - GetIconImage, 411
  - GetPixmap, 412
  - IconImageGenerator, 411
  - SetOutputDimensions, 412
  - SetOutsideValuePixel, 412
  - SetPixelMinMax, 412
  - SetPixmap, 412
- gdcmm::Image, 413
  - ~Image, 415
  - GetDirectionCosines, 415
  - GetIntercept, 415
  - GetOrigin, 415
  - GetSlope, 416
  - GetSpacing, 416
  - Image, 415
  - Print, 416
  - SetDirectionCosines, 416
  - SetIntercept, 416
  - SetOrigin, 416
  - SetSlope, 416
  - SetSpacing, 416
- gdcmm::ImageApplyLookupTable, 417
  - ~ImageApplyLookupTable, 419
  - Apply, 419
  - ImageApplyLookupTable, 419
- gdcmm::ImageChangePhotometricInterpretation, 419
  - ~ImageChangePhotometricInterpretation, 421
  - Change, 421
  - ChangeMonochrome, 421
  - GetPhotometricInterpretation, 421
  - ImageChangePhotometricInterpretation, 421
  - RGB2YBR, 421
  - SetPhotometricInterpretation, 421
  - YBR2RGB, 421
- gdcmm::ImageChangePlanarConfiguration, 422
  - ~ImageChangePlanarConfiguration, 424
  - Change, 424
  - GetPlanarConfiguration, 424
  - ImageChangePlanarConfiguration, 424
  - RGBPixelsToRGBPlanes, 424
  - RGBPlanesToRGBPixels, 424
  - SetPlanarConfiguration, 424
- gdcmm::ImageChangeTransferSyntax, 425
  - ~ImageChangeTransferSyntax, 427
  - Change, 427
  - GetTransferSyntax, 427
  - ImageChangeTransferSyntax, 427
  - SetCompressIconImage, 427
  - SetForce, 428
  - SetTransferSyntax, 428
  - SetUserCodec, 428
  - TryJPEG2000Codec, 428
  - TryJPEGCodec, 428
  - TryJPEGLSCodec, 428
  - TryRAWCodec, 428
  - TryRLECodec, 428
- gdcmm::ImageCodec, 429
  - ~ImageCodec, 431
  - AppendFrameEncode, 431
  - AppendRowEncode, 431
  - CanCode, 431
  - CanDecode, 432
  - Clone, 432
  - Decode, 432
  - DecodeByStreams, 432
  - Dimensions, 435
  - DoByteSwap, 432
  - DoInvertMonochrome, 432
  - DoOverlayCleanup, 432
  - DoPaddedCompositePixelCode, 432
  - DoPlanarConfiguration, 432
  - DoSimpleCopy, 432
  - DoYBR, 432
  - FileChangeTransferSyntax, 434
  - GetDimensions, 432
  - GetHeaderInfo, 433
  - GetLUT, 433
  - GetLossyFlag, 433
  - GetNeedByteSwap, 433
  - GetNumberOfDimensions, 433
  - GetPhotometricInterpretation, 433
  - GetPixelFormat, 433
  - GetPlanarConfiguration, 433
  - ImageChangePhotometricInterpretation, 434
  - ImageCodec, 431
  - IsFrameEncoder, 433
  - IsLossy, 433
  - IsRowEncoder, 433
  - IsValid, 433
  - LUT, 435
  - LUTPtr, 431
  - LossyFlag, 435
  - NeedByteSwap, 435
  - NeedOverlayCleanup, 435
  - NumberOfDimensions, 435
  - PF, 435
  - PI, 435
  - PlanarConfiguration, 435
  - RequestPaddedCompositePixelCode, 435
  - RequestPlanarConfiguration, 435

- SetDimensions, [433](#)
- SetLUT, [434](#)
- SetLossyFlag, [434](#)
- SetNeedByteSwap, [434](#)
- SetNeedOverlayCleanup, [434](#)
- SetNumberOfDimensions, [434](#)
- SetPhotometricInterpretation, [434](#)
- SetPixelFormat, [434](#)
- SetPlanarConfiguration, [434](#)
- StartEncode, [434](#)
- StopEncode, [434](#)
- gdcmm::ImageConverter, [435](#)
  - ~ImageConverter, [436](#)
  - Convert, [436](#)
  - GetOutput, [436](#)
  - ImageConverter, [436](#)
  - SetInput, [436](#)
- gdcmm::ImageFragmentSplitter, [436](#)
  - ~ImageFragmentSplitter, [439](#)
  - GetFragmentSizeMax, [439](#)
  - ImageFragmentSplitter, [439](#)
  - SetForce, [439](#)
  - SetFragmentSizeMax, [439](#)
  - Split, [439](#)
- gdcmm::ImageHelper, [439](#)
  - ComputeMediaStorageFromModality, [440](#)
  - ComputeSpacingFromImagePositionPatient, [440](#)
  - GetDimensionsValue, [441](#)
  - GetDirectionCosinesFromDataSet, [441](#)
  - GetDirectionCosinesValue, [441](#)
  - GetForcePixelSpacing, [441](#)
  - GetForceRescaleInterceptSlope, [441](#)
  - GetLUT, [441](#)
  - GetOriginValue, [441](#)
  - GetPhotometricInterpretationValue, [441](#)
  - GetPixelFormatValue, [441](#)
  - GetPlanarConfigurationValue, [441](#)
  - GetPointerFromElement, [441](#)
  - GetRescaleInterceptSlopeValue, [441](#)
  - GetSpacingTagFromMediaStorage, [442](#)
  - GetSpacingValue, [442](#)
  - GetZSpacingTagFromMediaStorage, [442](#)
  - SetDimensionsValue, [442](#)
  - SetDirectionCosinesValue, [442](#)
  - SetForcePixelSpacing, [442](#)
  - SetForceRescaleInterceptSlope, [442](#)
  - SetOriginValue, [442](#)
  - SetRescaleInterceptSlopeValue, [442](#)
  - SetSpacingValue, [442](#)
- gdcmm::ImageReader, [443](#)
  - ~ImageReader, [445](#)
  - GetImage, [445](#)
  - ImageReader, [445](#)
  - Read, [445](#)
- ReadACRNEMAIImage, [446](#)
- ReadImage, [446](#)
- gdcmm::ImageRegionReader, [446](#)
  - ~ImageRegionReader, [448](#)
  - ComputeBufferLength, [448](#)
  - GetRegion, [448](#)
  - ImageRegionReader, [448](#)
  - Read, [448](#)
  - ReadInformation, [448](#)
  - ReadIntoBuffer, [448](#)
  - SetRegion, [449](#)
- gdcmm::ImageToImageFilter, [449](#)
  - ~ImageToImageFilter, [451](#)
  - GetInput, [451](#)
  - GetOutput, [451](#)
  - ImageToImageFilter, [450](#)
- gdcmm::ImageWriter, [451](#)
  - ~ImageWriter, [453](#)
  - GetImage, [453](#)
  - ImageWriter, [453](#)
  - Write, [453](#)
- gdcmm::ImplicitDataElement, [456](#)
  - GetLength, [457](#)
  - Read, [457](#)
  - ReadPreValue, [457](#)
  - ReadValue, [457](#)
  - ReadValueWithLength, [457](#)
  - ReadWithLength, [457](#)
  - Write, [457](#)
- gdcmm::InitializeEvent, [458](#)
- gdcmm::Item, [467](#)
  - Clear, [469](#)
  - FindDataElement, [469](#)
  - GetDataElement, [469](#)
  - GetLength, [469](#)
  - GetNestedDataSet, [469](#)
  - InsertDataElement, [470](#)
  - Item, [469](#)
  - operator<<, [470](#)
  - Read, [470](#)
  - SetNestedDataSet, [470](#)
  - Write, [470](#)
- gdcmm::IterationEvent, [470](#)
- gdcmm::JPEG12Codec, [471](#)
  - ~JPEG12Codec, [473](#)
  - DecodeByStreams, [473](#)
  - EncodeBuffer, [473](#)
  - GetHeaderInfo, [473](#)
  - InternalCode, [473](#)
  - IsStateSuspension, [473](#)
  - JPEG12Codec, [473](#)
- gdcmm::JPEG16Codec, [474](#)
  - ~JPEG16Codec, [475](#)
  - DecodeByStreams, [475](#)

- EncodeBuffer, [475](#)
- GetHeaderInfo, [475](#)
- InternalCode, [475](#)
- IsStateSuspension, [475](#)
- JPEG16Codec, [475](#)
- gdcmm::JPEG2000Codec, [476](#)
  - ~JPEG2000Codec, [477](#)
  - AppendFrameEncode, [478](#)
  - AppendRowEncode, [478](#)
  - Bitmap, [479](#)
  - CanCode, [478](#)
  - CanDecode, [478](#)
  - Clone, [478](#)
  - Code, [478](#)
  - Decode, [478](#)
  - DecodeByStreams, [478](#)
  - DecodeExtent, [478](#)
  - GetHeaderInfo, [478](#)
  - GetQuality, [479](#)
  - GetRate, [479](#)
  - ImageRegionReader, [479](#)
  - IsFrameEncoder, [479](#)
  - IsRowEncoder, [479](#)
  - JPEG2000Codec, [477](#)
  - SetNumberOfResolutions, [479](#)
  - SetQuality, [479](#)
  - SetRate, [479](#)
  - SetReversible, [479](#)
  - SetTileSize, [479](#)
  - StartEncode, [479](#)
  - StopEncode, [479](#)
- gdcmm::JPEG8Codec, [479](#)
  - ~JPEG8Codec, [481](#)
  - DecodeByStreams, [481](#)
  - EncodeBuffer, [481](#)
  - GetHeaderInfo, [481](#)
  - InternalCode, [481](#)
  - IsStateSuspension, [481](#)
  - JPEG8Codec, [481](#)
- gdcmm::JPEGCodec, [482](#)
  - ~JPEGCodec, [484](#)
  - AppendFrameEncode, [484](#)
  - AppendRowEncode, [484](#)
  - BitSample, [486](#)
  - CanCode, [484](#)
  - CanDecode, [484](#)
  - Clone, [484](#)
  - Code, [484](#)
  - ComputeOffsetTable, [485](#)
  - Decode, [485](#)
  - DecodeByStreams, [485](#)
  - DecodeExtent, [485](#)
  - EncodeBuffer, [485](#)
  - GetHeaderInfo, [485](#)
  - GetLossless, [485](#)
  - GetQuality, [485](#)
  - ImageRegionReader, [486](#)
  - IsFrameEncoder, [485](#)
  - IsRowEncoder, [485](#)
  - IsStateSuspension, [485](#)
  - IsValid, [486](#)
  - JPEGCodec, [484](#)
  - Quality, [486](#)
  - SetBitSample, [486](#)
  - SetLossless, [486](#)
  - SetPixelFormat, [486](#)
  - SetQuality, [486](#)
  - StartEncode, [486](#)
  - StopEncode, [486](#)
- gdcmm::JPEGLSCodec, [486](#)
  - ~JPEGLSCodec, [488](#)
  - AppendFrameEncode, [488](#)
  - AppendRowEncode, [488](#)
  - CanCode, [489](#)
  - CanDecode, [489](#)
  - Clone, [489](#)
  - Code, [489](#)
  - Decode, [489](#)
  - DecodeExtent, [489](#)
  - GetBufferLength, [489](#)
  - GetHeaderInfo, [489](#)
  - GetLossless, [489](#)
  - ImageRegionReader, [490](#)
  - IsFrameEncoder, [489](#)
  - IsRowEncoder, [489](#)
  - JPEGLSCodec, [488](#)
  - SetBufferLength, [490](#)
  - SetLossless, [490](#)
  - SetLossyError, [490](#)
  - StartEncode, [490](#)
  - StopEncode, [490](#)
- gdcmm::JSON, [490](#)
  - ~JSON, [491](#)
  - Code, [491](#)
  - Decode, [491](#)
  - GetPrettyPrint, [491](#)
  - JSON, [491](#)
  - PrettyPrintOff, [491](#)
  - PrettyPrintOn, [491](#)
  - SetPrettyPrint, [491](#)
- gdcmm::KAKADUCodec, [491](#)
  - ~KAKADUCodec, [493](#)
  - CanCode, [493](#)
  - CanDecode, [493](#)
  - Clone, [493](#)
  - Code, [493](#)
  - Decode, [493](#)
  - KAKADUCodec, [493](#)

- gdcmm::LO, 493
  - const\_iterator, 495
  - const\_reference, 495
  - const\_reverse\_iterator, 495
  - difference\_type, 495
  - IsValid, 495
  - iterator, 495
  - LO, 495
  - pointer, 495
  - reference, 495
  - reverse\_iterator, 495
  - size\_type, 495
  - Superclass, 495
  - value\_type, 495
- gdcmm::LookupTable, 496
  - ~LookupTable, 498
  - Allocate, 498
  - BLUE, 498
  - BitSample, 500
  - Clear, 498
  - Decode, 498
  - GRAY, 498
  - GREEN, 498
  - GetBitSample, 499
  - GetBufferAsRGBA, 499
  - GetLUT, 499
  - GetLUTDescriptor, 499
  - GetLUTLength, 499
  - GetPointer, 499
  - IncompleteLUT, 500
  - InitializeBlueLUT, 499
  - InitializeGreenLUT, 499
  - InitializeLUT, 499
  - InitializeRedLUT, 499
  - Initialized, 499
  - Internal, 500
  - LookupTable, 498
  - LookupTableType, 498
  - Print, 499
  - RED, 498
  - SetBlueLUT, 500
  - SetGreenLUT, 500
  - SetLUT, 500
  - SetRedLUT, 500
  - UNKNOWN, 498
  - WriteBufferAsRGBA, 500
- gdcmm::MD5, 504
  - ~MD5, 505
  - Compute, 505
  - ComputeFile, 505
  - MD5, 505
- gdcmm::Macro, 501
  - AddMacroEntry, 502
  - ArrayIncludeMacrosType, 501
  - Clear, 502
  - FindMacroEntry, 502
  - GetMacroEntry, 502
  - GetName, 502
  - Macro, 501
  - MapModuleEntry, 501
  - operator<<, 502
  - SetName, 502
  - Verify, 502
- gdcmm::Macros, 502
  - AddMacro, 503
  - Clear, 503
  - GetMacro, 503
  - IsEmpty, 503
  - Macros, 503
  - ModuleMapType, 503
  - operator<<, 503
- gdcmm::MediaStorage, 505
  - AmbulatoryECGWaveformStorage, 509
  - Audio, 510
  - BasicTextSR, 509
  - BasicVoiceAudioWaveformStorage, 509
  - BreastTomosynthesisImageStorage, 510
  - CSANonImageStorage, 509
  - CTImageStorage, 508
  - CardiacElectrophysiologyWaveformStorage, 509
  - ComprehensiveSR, 509
  - ComputedRadiographyImageStorage, 508
  - DetachedPatientManagementSOPClass, 509
  - DetachedStudyManagementSOPClass, 509
  - DetachedVisitManagementSOPClass, 509
  - DigitalIntraoralXRayImageStorageForProcessing, 508
  - DigitalIntraoralXrayImageStorageForPresentation, 508
  - DigitalMammographyImageStorageForPresentation, 508
  - DigitalMammographyImageStorageForProcessing, 508
  - DigitalXRayImageStorageForPresentation, 508
  - DigitalXRayImageStorageForProcessing, 508
  - EncapsulatedCDASStorage, 509
  - EncapsulatedPDFStorage, 509
  - EnhancedCTImageStorage, 508
  - EnhancedMRIImageStorage, 508
  - EnhancedSR, 509
  - EnhancedUSVolumeStorage, 510
  - EnhancedXAImageStorage, 510
  - FujiPrivateCRLImageStorage, 510
  - GEPrivate3DModelStorage, 509
  - GeneralECGWaveformStorage, 509
  - GeneralElectricMagneticResonanceImageStorage, 509
  - GetMSString, 511

- GetMSType, 511
- GetModality, 511
- GetModalityDimension, 511
- GetNumberOfMSString, 511
- GetNumberOfMSType, 511
- GetNumberOfModality, 511
- GetString, 511
- GrayscaleSoftcopyPresentationStateStorageSOP↔  
Class, 509
- GuessFromModality, 511
- HangingProtocolStorage, 510
- HardcopyGrayscaleImageStorage, 509
- HemodynamicWaveformStorage, 509
- IsImage, 511
- IsUndefined, 512
- KeyObjectSelectionDocument, 509
- LeadECGWaveformStorage, 509
- MRImageStorage, 508
- MRSpectroscopyStorage, 508
- MS\_END, 510
- MSType, 508
- MammographyCADSR, 509
- MediaStorage, 511
- MediaStorageDirectoryStorage, 508
- ModalityPerformedProcedureStepSOPClass, 510
- MultiframeGrayscaleByteSecondaryCaptureImage↔  
Storage, 508
- MultiframeGrayscaleWordSecondaryCapture↔  
ImageStorage, 508
- MultiframeSingleBitSecondaryCaptureImage↔  
Storage, 508
- MultiframeTrueColorSecondaryCaptureImage↔  
Storage, 509
- NoObject, 510
- NuclearMedicineImageStorage, 509
- NuclearMedicineImageStorageRetired, 508
- ObjectEnd, 510
- ObjectType, 510
- operator MSType, 512
- operator<<, 512
- OphthalmicPhotography8BitImageStorage, 510
- OphthalmicTomographyImageStorage, 510
- PDF, 510
- PETImageStorage, 509
- Philips3D, 509
- PhilipsPrivateMRSyntheticImageStorage, 510
- RTDoseStorage, 509
- RTImageStorage, 509
- RTIonBeamsTreatmentRecordStorage, 510
- RTIonPlanStorage, 510
- RTPlanStorage, 509
- RTStructureSetStorage, 509
- RTTreatmentSummaryRecordStorage, 510
- RawDataStorage, 509
- SecondaryCaptureImageStorage, 508
- Segmentation, 510
- SegmentationStorage, 510
- SetFromDataSet, 512
- SetFromFile, 512
- SetFromHeader, 512
- SetFromModality, 512
- SetFromSourceImageSequence, 512
- SpacialFiducialsStorage, 509
- SpacialRegistrationStorage, 509
- StandaloneCurveStorage, 509
- StandaloneModalityLUTStorage, 509
- StandaloneOverlayStorage, 509
- StandaloneVOILUTStorage, 509
- StudyComponentManagementSOPClass, 509
- SurfaceSegmentationStorage, 510
- ToshibaPrivateDataStorage, 509
- URI, 510
- UltrasoundImageStorage, 508
- UltrasoundImageStorageRetired, 508
- UltrasoundMultiFrameImageStorage, 508
- UltrasoundMultiFrameImageStorageRetired, 508
- VLEndoscopicImageStorage, 510
- VLMicroscopicImageStorage, 510
- VLPhotographicImageStorage, 510
- VLWholeSlideMicroscopyImageStorage, 510
- Video, 510
- VideoEndoscopicImageStorage, 509
- Waveform, 510
- XRay3DAngiographicImageStorage, 510
- XRayAngiographicBiPlaneImageStorageRetired, 509
- XRayAngiographicImageStorage, 509
- XRayRadiationDoseSR, 510
- XRayRadiofluoroscopicImageStorage, 509
- gdcmmemberCommand
  - ~MemberCommand, 515
  - Execute, 515
  - m\_ConstMemberFunction, 515
  - m\_MemberFunction, 516
  - m\_This, 516
  - MemberCommand, 515
  - New, 515
  - Self, 514
  - SetCallbackFunction, 515
  - TConstMemberFunctionPointer, 514
  - TMemberFunctionPointer, 514
- gdcmmemberCommand< T >, 512
- gdcmmeshPrimitive, 516
  - ~MeshPrimitive, 519
  - AddPrimitiveData, 519
  - EDGE, 518
  - FACET, 518
  - GetMPType, 519
  - GetMPTypeString, 519



- GetNumberOfPrimitivesData, [519](#)
- GetPrimitiveData, [519](#)
- GetPrimitiveType, [519](#)
- GetPrimitivesData, [519](#)
- LINE, [518](#)
- MPTType, [518](#)
- MPTType\_END, [518](#)
- MeshPrimitive, [519](#)
- PrimitiveData, [519](#)
- PrimitiveType, [519](#)
- PrimitivesData, [518](#)
- SetPrimitiveData, [519](#)
- SetPrimitiveType, [519](#)
- SetPrimitivesData, [519](#)
- TRIANGLE, [518](#)
- TRIANGLE\_FAN, [518](#)
- TRIANGLE\_STRIP, [518](#)
- VERTEX, [518](#)
- gdcmm::ModifiedEvent, [519](#)
- gdcmm::Module, [521](#)
  - AddMacro, [522](#)
  - AddModuleEntry, [522](#)
  - ArrayIncludeMacrosType, [521](#)
  - Clear, [522](#)
  - FindModuleEntryInMacros, [522](#)
  - GetModuleEntryInMacros, [522](#)
  - GetName, [522](#)
  - MapModuleEntry, [521](#)
  - Module, [522](#)
  - operator<<, [522](#)
  - SetName, [522](#)
  - Verify, [522](#)
- gdcmm::ModuleEntry, [522](#)
  - ~ModuleEntry, [524](#)
  - DataElementType, [525](#)
  - Description, [524](#)
  - DescriptionField, [525](#)
  - GetDescription, [524](#)
  - GetName, [524](#)
  - GetType, [524](#)
  - ModuleEntry, [524](#)
  - Name, [525](#)
  - operator<<, [525](#)
  - SetDescription, [525](#)
  - SetName, [525](#)
  - SetType, [525](#)
- gdcmm::Modules, [525](#)
  - AddModule, [526](#)
  - Clear, [526](#)
  - GetModule, [526](#)
  - IsEmpty, [526](#)
  - ModuleMapType, [526](#)
  - Modules, [526](#)
  - operator<<, [526](#)
- gdcmm::MovePatientRootQuery, [527](#)
  - GetAbstractSyntaxUID, [528](#)
  - GetTagListByLevel, [528](#)
  - InitializeDataSet, [528](#)
  - MovePatientRootQuery, [528](#)
  - QueryFactory, [529](#)
  - ValidateQuery, [528](#)
- gdcmm::MoveStudyRootQuery, [529](#)
  - GetAbstractSyntaxUID, [530](#)
  - GetTagListByLevel, [530](#)
  - InitializeDataSet, [530](#)
  - MoveStudyRootQuery, [530](#)
  - QueryFactory, [531](#)
  - ValidateQuery, [530](#)
- gdcmm::NestedModuleEntries, [531](#)
  - AddModuleEntry, [533](#)
  - GetModuleEntry, [533](#)
  - GetNumberOfModuleEntries, [533](#)
  - NestedModuleEntries, [533](#)
  - operator<<, [533](#)
  - SizeType, [533](#)
- gdcmm::NoEvent, [533](#)
- gdcmm::Object, [534](#)
  - ~Object, [536](#)
  - Object, [536](#)
  - operator<<, [536](#)
  - operator=, [536](#)
  - Print, [536](#)
  - Register, [536](#)
  - SmartPointer, [536](#)
  - UnRegister, [536](#)
- gdcmm::OpenSSLCryptoFactory, [537](#)
  - CreateCMSProvider, [538](#)
  - InitOpenSSL, [538](#)
  - OpenSSLCryptoFactory, [538](#)
- gdcmm::OpenSSLCryptographicMessageSyntax, [538](#)
  - ~OpenSSLCryptographicMessageSyntax, [539](#)
  - Decrypt, [539](#)
  - Encrypt, [539](#)
  - GetCipherType, [540](#)
  - OpenSSLCryptographicMessageSyntax, [539](#)
  - ParseCertificateFile, [540](#)
  - ParseKeyFile, [540](#)
  - SetCipherType, [540](#)
  - SetPassword, [540](#)
- gdcmm::OpenSSLP7CryptoFactory, [540](#)
  - CreateCMSProvider, [541](#)
  - OpenSSLP7CryptoFactory, [541](#)
- gdcmm::OpenSSLP7CryptographicMessageSyntax, [542](#)
  - ~OpenSSLP7CryptographicMessageSyntax, [543](#)
  - Decrypt, [543](#)
  - Encrypt, [543](#)
  - GetCipherType, [543](#)
  - OpenSSLP7CryptographicMessageSyntax, [543](#)



- ParseCertificateFile, [543](#)
- ParseKeyFile, [544](#)
- SetCipherType, [544](#)
- SetPassword, [544](#)
- gdcmm::Orientation, [544](#)
  - ~Orientation, [545](#)
  - AXIAL, [545](#)
  - CORONAL, [545](#)
  - GetLabel, [545](#)
  - GetMajorAxisFromPatientRelativeDirectionCosine, [545](#)
  - GetObliquityThresholdCosineValue, [546](#)
  - GetType, [546](#)
  - OBLIQUE, [545](#)
  - operator<<, [546](#)
  - Orientation, [545](#)
  - OrientationType, [545](#)
  - Print, [546](#)
  - SAGITTAL, [545](#)
  - SetObliquityThresholdCosineValue, [546](#)
  - UNKNOWN, [545](#)
- gdcmm::Overlay, [546](#)
  - ~Overlay, [549](#)
  - Decompress, [549](#)
  - GetBitPosition, [549](#)
  - GetBitsAllocated, [549](#)
  - GetColumns, [549](#)
  - GetDescription, [549](#)
  - GetGroup, [550](#)
  - GetOrigin, [550](#)
  - GetOverlayData, [550](#)
  - GetOverlayTypeAsString, [550](#)
  - GetOverlayTypeFromString, [550](#)
  - GetRows, [550](#)
  - GetType, [550](#)
  - GetTypeAsEnum, [550](#)
  - GetUnpackBuffer, [550](#)
  - GetUnpackBufferLength, [550](#)
  - GrabOverlayFromPixelData, [550](#)
  - Graphics, [549](#)
  - Invalid, [549](#)
  - IsEmpty, [550](#)
  - IsInPixelData, [550](#)
  - IsZero, [551](#)
  - Overlay, [549](#)
  - OverlayType, [549](#)
  - Print, [551](#)
  - ROI, [549](#)
  - SetBitPosition, [551](#)
  - SetBitsAllocated, [551](#)
  - SetColumns, [551](#)
  - SetDescription, [551](#)
  - SetFrameOrigin, [551](#)
  - SetGroup, [551](#)
  - SetNumberOfFrames, [551](#)
  - SetOrigin, [551](#)
  - SetOverlay, [551](#)
  - SetRows, [552](#)
  - SetType, [552](#)
  - Update, [552](#)
- gdcmm::PDBElement, [559](#)
  - GetName, [560](#)
  - GetValue, [560](#)
  - NameField, [560](#)
  - operator<<, [560](#)
  - operator==, [560](#)
  - PDBElement, [560](#)
  - SetName, [560](#)
  - SetValue, [560](#)
  - ValueField, [560](#)
- gdcmm::PDBHeader, [561](#)
  - ~PDBHeader, [562](#)
  - FindPDBElementByName, [562](#)
  - GetPDBEEnd, [562](#)
  - GetPDBElementByName, [562](#)
  - GetPDBInfoTag, [562](#)
  - LoadFromDataElement, [562](#)
  - operator<<, [562](#)
  - PDBHeader, [562](#)
  - Print, [562](#)
- gdcmm::PDFCodec, [562](#)
  - ~PDFCodec, [564](#)
  - CanCode, [564](#)
  - CanDecode, [564](#)
  - Decode, [564](#)
  - PDFCodec, [564](#)
- gdcmm::PGXCodec, [567](#)
  - ~PGXCodec, [568](#)
  - CanCode, [568](#)
  - CanDecode, [568](#)
  - Clone, [568](#)
  - GetHeaderInfo, [569](#)
  - PGXCodec, [568](#)
  - Read, [569](#)
  - Write, [569](#)
- gdcmm::PNMCodec, [588](#)
  - ~PNMCodec, [590](#)
  - CanCode, [590](#)
  - CanDecode, [590](#)
  - Clone, [590](#)
  - GetBufferLength, [590](#)
  - GetHeaderInfo, [590](#)
  - PNMCodec, [590](#)
  - Read, [590](#)
  - SetBufferLength, [590](#)
  - Write, [590](#)
- gdcmm::PVRGCodec, [610](#)
  - ~PVRGCodec, [611](#)

- CanCode, [611](#)
- CanDecode, [611](#)
- Clone, [611](#)
- Code, [612](#)
- Decode, [612](#)
- PVRGCodec, [611](#)
- SetLossyFlag, [612](#)
- gdcmm::ParseException, [552](#)
  - ~ParseException, [553](#)
  - GetLastElement, [553](#)
  - operator=, [553](#)
  - ParseException, [553](#)
  - SetLastElement, [553](#)
- gdcmm::Parser, [554](#)
  - ~Parser, [555](#)
  - DuplicateAttributeError, [555](#)
  - EndElementHandler, [555](#)
  - ErrorType, [555](#)
  - GetBuffer, [555](#)
  - GetCurrentByteIndex, [555](#)
  - GetErrorCode, [555](#)
  - GetErrorString, [555](#)
  - GetUserData, [555](#)
  - JunkAfterDocElementError, [555](#)
  - NoElementsError, [555](#)
  - NoError, [555](#)
  - NoMemoryError, [555](#)
  - Parse, [555](#)
  - ParseBuffer, [556](#)
  - Parser, [555](#)
  - Process, [556](#)
  - SetElementHandler, [556](#)
  - SetUserData, [556](#)
  - StartElementHandler, [555](#)
  - SyntaxError, [555](#)
  - TagMismatchError, [555](#)
  - UndefinedEntityError, [555](#)
  - UnexpectedStateError, [555](#)
- gdcmm::Patient, [556](#)
  - Patient, [556](#)
- gdcmm::PersonName, [566](#)
  - Component, [566](#)
  - GetMaxLength, [566](#)
  - GetNumberOfComponents, [566](#)
  - MaxLength, [567](#)
  - MaxNumberOfComponents, [567](#)
  - Padding, [567](#)
  - Print, [566](#)
  - Separator, [567](#)
  - SetBlob, [566](#)
  - SetComponents, [566](#)
- gdcmm::PhotometricInterpretation, [569](#)
  - ARGB, [570](#)
  - CMYK, [570](#)
  - GetPIString, [571](#)
  - GetPIType, [571](#)
  - GetSamplesPerPixel, [571](#)
  - GetString, [571](#)
  - GetType, [571](#)
  - HSV, [570](#)
  - IsLossless, [571](#)
  - IsLossy, [571](#)
  - IsRetired, [571](#)
  - IsSameColorSpace, [571](#)
  - MONOCHROME1, [570](#)
  - MONOCHROME2, [570](#)
  - operator PIType, [571](#)
  - operator<<, [571](#)
  - PALETTE\_COLOR, [570](#)
  - PI\_END, [570](#)
  - PIType, [570](#)
  - PhotometricInterpretation, [571](#)
  - RGB, [570](#)
  - UNKNOWN, [570](#)
  - YBR\_FULL, [570](#)
  - YBR\_FULL\_422, [570](#)
  - YBR\_ICT, [570](#)
  - YBR\_PARTIAL\_420, [570](#)
  - YBR\_PARTIAL\_422, [570](#)
  - YBR\_RCT, [570](#)
- gdcmm::PixelFormat, [571](#)
  - Bitmap, [576](#)
  - FLOAT16, [574](#)
  - FLOAT32, [574](#)
  - FLOAT64, [574](#)
  - GetBitsAllocated, [574](#)
  - GetBitsStored, [574](#)
  - GetHighBit, [574](#)
  - GetMax, [574](#)
  - GetMin, [574](#)
  - GetPixelRepresentation, [575](#)
  - GetPixelSize, [575](#)
  - GetSamplesPerPixel, [575](#)
  - GetScalarType, [575](#)
  - GetScalarTypeAsString, [575](#)
  - INT12, [574](#)
  - INT16, [574](#)
  - INT32, [574](#)
  - INT64, [574](#)
  - INT8, [573](#)
  - IsCompatible, [575](#)
  - IsValid, [575](#)
  - operator ScalarType, [575](#)
  - operator!=, [575](#)
  - operator<<, [576](#)
  - operator==, [575](#), [576](#)
  - PixelFormat, [574](#)
  - Print, [576](#)

- SINGLEBIT, [574](#)
- ScalarType, [573](#)
- SetBitsAllocated, [576](#)
- SetBitsStored, [576](#)
- SetHighBit, [576](#)
- SetPixelRepresentation, [576](#)
- SetSamplesPerPixel, [576](#)
- SetScalarType, [576](#)
- UINT12, [573](#)
- UINT16, [574](#)
- UINT32, [574](#)
- UINT64, [574](#)
- UINT8, [573](#)
- UNKNOWN, [574](#)
- Validate, [576](#)
- gdcm::Pixmap, [577](#)
  - ~Pixmap, [578](#)
  - AreOverlaysInPixelData, [578](#)
  - Curves, [579](#)
  - GetCurve, [578](#), [579](#)
  - GetIconImage, [579](#)
  - GetNumberOfCurves, [579](#)
  - GetNumberOfOverlays, [579](#)
  - GetOverlay, [579](#)
  - Icon, [579](#)
  - Overlays, [579](#)
  - Pixmap, [578](#)
  - Print, [579](#)
  - RemoveOverlay, [579](#)
  - SetIconImage, [579](#)
  - SetNumberOfCurves, [579](#)
  - SetNumberOfOverlays, [579](#)
- gdcm::PixmapReader, [579](#)
  - ~PixmapReader, [582](#)
  - GetPixmap, [582](#)
  - PixelData, [583](#)
  - PixmapReader, [582](#)
  - Read, [582](#)
  - ReadACRNEMAImage, [582](#)
  - ReadImage, [582](#)
  - ReadImageInternal, [582](#)
- gdcm::PixmapToPixmapFilter, [583](#)
  - ~PixmapToPixmapFilter, [584](#)
  - GetInput, [585](#)
  - GetOutput, [585](#)
  - GetOutputAsPixmap, [585](#)
  - PixmapToPixmapFilter, [584](#)
- gdcm::PixmapWriter, [585](#)
  - ~PixmapWriter, [587](#)
  - DolconImage, [587](#)
  - GetImage, [587](#)
  - GetPixmap, [587](#)
  - PixelData, [588](#)
  - PixmapWriter, [587](#)
  - PrepareWrite, [587](#)
  - SetImage, [587](#)
  - SetPixmap, [588](#)
  - Write, [588](#)
- gdcm::Preamble, [591](#)
  - ~Preamble, [592](#)
  - Clear, [592](#)
  - Create, [592](#)
  - GetInternal, [592](#)
  - GetLength, [592](#)
  - IsEmpty, [592](#)
  - IsValid, [592](#)
  - operator<<, [592](#)
  - operator=, [592](#)
  - Preamble, [592](#)
  - Print, [592](#)
  - Read, [592](#)
  - Remove, [592](#)
  - Valid, [592](#)
  - Write, [592](#)
- gdcm::PresentationContext, [592](#)
  - AddTransferSyntax, [593](#)
  - GetAbstractSyntax, [593](#)
  - GetNumberOfTransferSyntaxes, [593](#)
  - GetPresentationContextID, [593](#)
  - GetTransferSyntax, [594](#)
  - operator==, [594](#)
  - PresentationContext, [593](#)
  - Print, [594](#)
  - SetAbstractSyntax, [594](#)
  - SetPresentationContextID, [594](#)
  - SizeType, [593](#)
  - TransferSyntaxArrayType, [593](#)
- gdcm::PresentationContextGenerator, [595](#)
  - AddPresentationContext, [596](#)
  - GenerateFromFilenames, [596](#)
  - GenerateFromUID, [597](#)
  - GetDefaultTransferSyntax, [597](#)
  - GetPresentationContexts, [597](#)
  - PresentationContextArrayType, [596](#)
  - PresentationContextGenerator, [596](#)
  - SetDefaultTransferSyntax, [597](#)
  - SetMergeModeToAbstractSyntax, [597](#)
  - SetMergeModeToTransferSyntax, [597](#)
  - SizeType, [596](#)
- gdcm::Printer, [601](#)
  - ~Printer, [603](#)
  - CONDENSED\_STYLE, [603](#)
  - F, [604](#)
  - GetPrintStyle, [603](#)
  - MaxPrintLength, [604](#)
  - Print, [603](#)
  - PrintDataElement, [603](#)
  - PrintDataSet, [603](#)

- PrintSQ, [603](#)
- PrintStyle, [604](#)
- PrintStyles, [603](#)
- Printer, [603](#)
- SetColor, [604](#)
- SetFile, [604](#)
- SetStyle, [604](#)
- VERBOSE\_STYLE, [603](#)
- XML, [603](#)
- gdcmm::PrivateDict, [604](#)
  - ~PrivateDict, [605](#)
  - AddDictEntry, [605](#)
  - Dicts, [605](#)
  - FindDictEntry, [605](#)
  - GetDictEntry, [605](#)
  - IsEmpty, [605](#)
  - LoadDefault, [605](#)
  - operator<<, [605](#)
  - PrintXML, [605](#)
  - PrivateDict, [605](#)
  - RemoveDictEntry, [605](#)
- gdcmm::PrivateTag, [606](#)
  - GetAsDataElement, [607](#)
  - GetOwner, [607](#)
  - operator<, [607](#)
  - operator<<, [607](#)
  - PrivateTag, [607](#)
  - ReadFromCommaSeparatedString, [607](#)
  - SetOwner, [607](#)
- gdcmm::ProgressEvent, [608](#)
  - ~ProgressEvent, [609](#)
  - CheckEvent, [609](#)
  - GetEventName, [609](#)
  - GetProgress, [609](#)
  - MakeObject, [609](#)
  - ProgressEvent, [609](#)
  - Self, [609](#)
  - SetProgress, [609](#)
  - Superclass, [609](#)
- gdcmm::PythonFilter, [612](#)
  - ~PythonFilter, [613](#)
  - GetFile, [613](#)
  - PythonFilter, [613](#)
  - SetDicts, [613](#)
  - SetFile, [613](#)
  - ToPyObject, [613](#)
  - UseDictAlways, [613](#)
- gdcmm::QueryBase, [613](#)
  - ~QueryBase, [614](#)
  - GetAllRequiredTags, [614](#)
  - GetAllTags, [614](#)
  - GetHierarchicalSearchTags, [614](#)
  - GetName, [614](#)
  - GetOptionalTags, [615](#)
  - GetQueryLevel, [615](#)
  - GetRequiredTags, [615](#)
  - GetUniqueTags, [615](#)
- gdcmm::QueryFactory, [615](#)
  - GetCharacterFromCurrentLocale, [616](#)
  - ListCharSets, [616](#)
  - ProduceCharacterSetDataElement, [616](#)
  - ProduceQuery, [616](#)
- gdcmm::QueryImage, [616](#)
  - GetHierarchicalSearchTags, [617](#)
  - GetName, [618](#)
  - GetOptionalTags, [618](#)
  - GetQueryLevel, [618](#)
  - GetRequiredTags, [618](#)
  - GetUniqueTags, [618](#)
- gdcmm::QueryPatient, [618](#)
  - GetHierarchicalSearchTags, [619](#)
  - GetName, [620](#)
  - GetOptionalTags, [620](#)
  - GetQueryLevel, [620](#)
  - GetRequiredTags, [620](#)
  - GetUniqueTags, [620](#)
- gdcmm::QuerySeries, [620](#)
  - GetHierarchicalSearchTags, [621](#)
  - GetName, [622](#)
  - GetOptionalTags, [622](#)
  - GetQueryLevel, [622](#)
  - GetRequiredTags, [622](#)
  - GetUniqueTags, [622](#)
- gdcmm::QueryStudy, [622](#)
  - GetHierarchicalSearchTags, [623](#)
  - GetName, [624](#)
  - GetOptionalTags, [624](#)
  - GetQueryLevel, [624](#)
  - GetRequiredTags, [624](#)
  - GetUniqueTags, [624](#)
- gdcmm::RAWCodec, [624](#)
  - ~RAWCodec, [626](#)
  - CanCode, [626](#)
  - CanDecode, [626](#)
  - Clone, [626](#)
  - Code, [626](#)
  - Decode, [626](#)
  - DecodeByStreams, [626](#)
  - DecodeBytes, [627](#)
  - GetHeaderInfo, [627](#)
  - RAWCodec, [626](#)
- gdcmm::RLECodec, [637](#)
  - ~RLECodec, [638](#)
  - AppendFrameEncode, [638](#)
  - AppendRowEncode, [639](#)
  - CanCode, [639](#)
  - CanDecode, [639](#)
  - Clone, [639](#)

- Code, [639](#)
- Decode, [639](#)
- DecodeByStreams, [639](#)
- DecodeExtent, [639](#)
- GetBufferLength, [639](#)
- GetHeaderInfo, [639](#)
- ImageRegionReader, [640](#)
- IsFrameEncoder, [639](#)
- IsRowEncoder, [640](#)
- RLECodec, [638](#)
- SetBufferLength, [640](#)
- SetLength, [640](#)
- StartEncode, [640](#)
- StopEncode, [640](#)
- gdcm::Reader, [627](#)
  - ~Reader, [630](#)
  - CanRead, [630](#)
  - F, [632](#)
  - GetFile, [630](#)
  - GetStreamCurrentPosition, [630](#)
  - GetStreamPtr, [630](#)
  - Read, [630](#)
  - ReadDataSet, [631](#)
  - ReadMetaInformation, [631](#)
  - ReadPreamble, [631](#)
  - ReadSelectedPrivateTags, [631](#)
  - ReadSelectedTags, [631](#)
  - ReadUpToTag, [631](#)
  - Reader, [630](#)
  - SetFile, [631](#)
  - SetFileName, [631](#)
  - SetStream, [632](#)
  - StreamImageReader, [632](#)
- gdcm::Region, [632](#)
  - ~Region, [633](#)
  - Area, [633](#)
  - Clone, [633](#)
  - ComputeBoundingBox, [633](#)
  - Empty, [633](#)
  - IsValid, [633](#)
  - Print, [634](#)
  - Region, [633](#)
- gdcm::Rescaler, [634](#)
  - ~Rescaler, [635](#)
  - ComputeInterceptSlopePixelType, [635](#)
  - ComputePixelTypeFromMinMax, [635](#)
  - GetIntercept, [636](#)
  - GetSlope, [636](#)
  - InverseRescale, [636](#)
  - InverseRescaleFunctionIntoBestFit, [636](#)
  - Rescale, [636](#)
  - RescaleFunctionIntoBestFit, [636](#)
  - Rescaler, [635](#)
  - SetIntercept, [636](#)
  - SetMinMaxForPixelType, [636](#)
  - SetPixelFormat, [636](#)
  - SetSlope, [636](#)
  - SetTargetPixelType, [636](#)
  - SetUseTargetPixelType, [636](#)
- gdcm::SHA1, [679](#)
  - ~SHA1, [680](#)
  - Compute, [680](#)
  - ComputeFile, [680](#)
  - SHA1, [680](#)
- gdcm::SOPClassUIDToIOD, [689](#)
  - const, [690](#)
  - GetIOD, [690](#)
  - GetIODFromSOPClassUID, [690](#)
  - GetNumberOfSOPClassToIOD, [690](#)
  - GetSOPClassUIDFromIOD, [690](#)
  - GetSOPClassUIDToIOD, [690](#)
  - GetSOPClassUIDToIODs, [690](#)
- gdcm::STATIC\_ASSERTION\_FAILURE< true >, [699](#)
  - value, [699](#)
- gdcm::STATIC\_ASSERTION\_FAILURE< x >, [699](#)
- gdcm::Scanner, [642](#)
  - ~Scanner, [645](#)
  - AddPrivateTag, [645](#)
  - AddSkipTag, [645](#)
  - AddTag, [645](#)
  - Begin, [645](#)
  - ClearSkipTags, [646](#)
  - ClearTags, [646](#)
  - ConstIterator, [645](#)
  - End, [646](#)
  - GetAllFilenamesFromTagToValue, [646](#)
  - GetFilenameFromTagToValue, [646](#)
  - GetFilenames, [646](#)
  - GetKeys, [646](#)
  - GetMapping, [646](#)
  - GetMappingFromTagToValue, [646](#)
  - GetMappings, [646](#)
  - GetOrderedValues, [646](#)
  - GetValue, [646](#)
  - GetValues, [647](#)
  - IsKey, [647](#)
  - MappingType, [645](#)
  - New, [647](#)
  - operator<<, [648](#)
  - Print, [647](#)
  - ProcessPublicTag, [647](#)
  - Scan, [647](#)
  - Scanner, [645](#)
  - TagToValue, [645](#)
  - TagToValueValueType, [645](#)
  - ValueType, [645](#)
- gdcm::Scanner::Itstr, [500](#)
  - operator(), [500](#)

- gdcmm::Segment, 648
  - ~Segment, 650
  - ALGOType, 650
  - ALGOType\_END, 650
  - AUTOMATIC, 650
  - AddSurface, 650
  - AnatomicRegion, 651
  - GetALGOType, 650
  - GetALGOTypeString, 650
  - GetAnatomicRegion, 650
  - GetPropertyCategory, 650
  - GetPropertyType, 650
  - GetSegmentAlgorithmName, 650
  - GetSegmentAlgorithmType, 651
  - GetSegmentDescription, 651
  - GetSegmentLabel, 651
  - GetSegmentNumber, 651
  - GetSurface, 651
  - GetSurfaceCount, 651
  - GetSurfaces, 651
  - MANUAL, 650
  - PropertyCategory, 651
  - PropertyType, 651
  - Segment, 650
  - SegmentAlgorithmName, 651
  - SegmentAlgorithmType, 651
  - SegmentDescription, 651
  - SegmentLabel, 651
  - SegmentNumber, 652
  - SetAnatomicRegion, 651
  - SetPropertyCategory, 651
  - SetPropertyType, 651
  - SetSegmentAlgorithmName, 651
  - SetSegmentAlgorithmType, 651
  - SetSegmentDescription, 651
  - SetSegmentLabel, 651
  - SetSegmentNumber, 651
  - SetSurfaceCount, 651
  - SurfaceCount, 652
  - SurfaceVector, 650
  - Surfaces, 652
- gdcmm::SegmentHelper, 136
- gdcmm::SegmentHelper::BasicCodedEntry, 204
  - BasicCodedEntry, 206
  - CM, 206
  - CSD, 206
  - CSV, 206
  - CV, 206
  - IsEmpty, 206
- gdcmm::SegmentReader, 654
  - ~SegmentReader, 656
  - GetSegments, 656
  - Read, 656
  - ReadSegment, 656
  - ReadSegments, 656
  - SegmentMap, 656
  - SegmentReader, 656
  - SegmentVector, 656
  - Segments, 656
- gdcmm::SegmentWriter, 656
  - ~SegmentWriter, 658
  - AddSegment, 658
  - GetNumberOfSegments, 658
  - GetSegment, 658
  - GetSegments, 658
  - PrepareWrite, 658
  - SegmentVector, 658
  - SegmentWriter, 658
  - Segments, 659
  - SetNumberOfSegments, 658
  - SetSegments, 658
  - Write, 658
- gdcmm::SegmentedPaletteColorLookupTable, 652
  - ~SegmentedPaletteColorLookupTable, 653
  - Print, 653
  - SegmentedPaletteColorLookupTable, 653
  - SetLUT, 653
- gdcmm::SequenceOfFragments, 659
  - AddFragment, 662
  - Begin, 662
  - Clear, 662
  - ComputeByteLength, 662
  - ComputeLength, 662
  - ConstIterator, 661
  - End, 662
  - FragmentVector, 661
  - GetBuffer, 662
  - GetFragBuffer, 662
  - GetFragment, 662
  - GetLength, 662
  - GetNumberOfFragments, 662
  - GetTable, 662, 663
  - Iterator, 661
  - New, 663
  - operator==, 663
  - Print, 663
  - Read, 663
  - ReadPreValue, 663
  - ReadValue, 663
  - SequenceOfFragments, 661
  - SetLength, 663
  - SizeType, 661
  - Write, 663
  - WriteBuffer, 663
- gdcmm::SequenceOfItems, 664
  - AddItem, 667
  - Begin, 667
  - Clear, 667

- ComputeLength, [667](#)
- ConstIterator, [667](#)
- End, [667](#)
- FindDataElement, [668](#)
- GetItem, [668](#)
- GetLength, [668](#)
- GetNumberOfItems, [668](#)
- IsUndefinedLength, [668](#)
- ItemVector, [667](#)
- Items, [669](#)
- Iterator, [667](#)
- New, [668](#)
- operator=, [668](#)
- operator==, [668](#)
- Print, [668](#)
- Read, [669](#)
- RemoveItemByIndex, [669](#)
- SequenceLengthField, [669](#)
- SequenceOfItems, [667](#)
- SetLength, [669](#)
- SetLengthToUndefined, [669](#)
- SetNumberOfItems, [669](#)
- SizeType, [667](#)
- Write, [669](#)
- gdcmm::SerieHelper, [670](#)
  - ~SerieHelper, [671](#)
  - AddFile, [672](#)
  - AddFileName, [672](#)
  - AddRestriction, [672](#)
  - Clear, [672](#)
  - CreateDefaultUniqueSeriesIdentifier, [672](#)
  - CreateUniqueSeriesIdentifier, [672](#)
  - FileNameOrdering, [672](#)
  - GetFirstSingleSerieUIDFileSet, [672](#)
  - GetNextSingleSerieUIDFileSet, [672](#)
  - ImagePositionPatientOrdering, [672](#)
  - ItFileSetHt, [672](#)
  - OrderFileList, [672](#)
  - SerieHelper, [671](#)
  - SerieRestrictions, [671](#)
  - SetDirectory, [672](#)
  - SetLoadMode, [672](#)
  - SetUseSeriesDetails, [672](#)
  - SingleSerieUIDFileSetHT, [672](#)
  - SingleSerieUIDFileSetmap, [671](#)
  - UserOrdering, [672](#)
- gdcmm::SerieHelper::Rule, [641](#)
  - elem, [642](#)
  - group, [642](#)
  - op, [642](#)
  - value, [642](#)
- gdcmm::Series, [672](#)
  - Series, [673](#)
- gdcmm::ServiceClassUser, [674](#)
  - ~ServiceClassUser, [676](#)
  - GetAETitle, [676](#)
  - GetCalledAETitle, [676](#)
  - GetTimeout, [677](#)
  - InitializeConnection, [677](#)
  - IsPresentationContextAccepted, [677](#)
  - New, [677](#)
  - SendEcho, [677](#)
  - SendFind, [677](#)
  - SendMove, [677](#)
  - SendStore, [677](#), [678](#)
  - ServiceClassUser, [676](#)
  - SetAETitle, [678](#)
  - SetCalledAETitle, [678](#)
  - SetHostname, [678](#)
  - SetPort, [678](#)
  - SetPortSCP, [678](#)
  - SetPresentationContexts, [678](#)
  - SetTimeout, [678](#)
  - StartAssociation, [679](#)
  - StopAssociation, [679](#)
- gdcmm::SimpleMemberCommand
  - ~SimpleMemberCommand, [683](#)
  - Execute, [683](#)
  - m\_MemberFunction, [683](#)
  - m\_This, [683](#)
  - New, [683](#)
  - Self, [682](#)
  - SetCallbackFunction, [683](#)
  - SimpleMemberCommand, [683](#)
  - TMemberFunctionPointer, [682](#)
- gdcmm::SimpleMemberCommand< T >, [680](#)
- gdcmm::SimpleSubjectWatcher, [684](#)
  - ~SimpleSubjectWatcher, [684](#)
  - EndFilter, [684](#)
  - ShowAbort, [685](#)
  - ShowAnonymization, [685](#)
  - ShowData, [685](#)
  - ShowDataSet, [685](#)
  - ShowFileName, [685](#)
  - ShowIteration, [685](#)
  - ShowProgress, [685](#)
  - SimpleSubjectWatcher, [684](#)
  - StartFilter, [685](#)
  - TestAbortOff, [685](#)
  - TestAbortOn, [685](#)
- gdcmm::SmartPointer
  - ~SmartPointer, [687](#)
  - GetPointer, [687](#)
  - operator ObjectType \*, [687](#)
  - operator\*, [688](#)
  - operator->, [688](#)
  - operator=, [688](#)
  - SmartPointer, [687](#)



- gdcmm::SmartPointer< ObjectType >, 685
- gdcmm::Sorter, 690
  - ~Sorter, 693
  - AddSelect, 693
  - FileNames, 694
  - GetFileNames, 693
  - operator<<, 694
  - Print, 693
  - Selection, 694
  - SelectionMap, 692
  - SetSortFunction, 693
  - Sort, 693
  - SortFunc, 694
  - SortFunction, 692
  - Sorter, 693
  - StableSort, 693
- gdcmm::Spacing, 694
  - ~Spacing, 695
  - CALIBRATED, 695
  - ComputePixelAspectRatioFromPixelSpacing, 695
  - DETECTOR, 695
  - MAGNIFIED, 695
  - Spacing, 695
  - SpacingType, 695
  - UNKNOWN, 695
- gdcmm::Spectroscopy, 696
  - Spectroscopy, 696
- gdcmm::SplitMosaicFilter, 696
  - ~SplitMosaicFilter, 697
  - ComputeMOSAICDimensions, 697
  - GetFile, 697
  - GetImage, 697
  - SetFile, 697
  - SetImage, 697
  - Split, 697
  - SplitMosaicFilter, 697
- gdcmm::StartEvent, 697
- gdcmm::StreamImageReader, 699
  - ~StreamImageReader, 700
  - CanReadImage, 700
  - DefinePixelExtent, 700
  - DefineProperBufferLength, 701
  - GetDimensionsValueForResolution, 701
  - GetFile, 701
  - Read, 701
  - ReadImageInformation, 701
  - SetFileName, 701
  - SetStream, 702
  - StreamImageReader, 700
- gdcmm::StreamImageWriter, 702
  - ~StreamImageWriter, 704
  - CanWriteFile, 704
  - DefinePixelExtent, 705
  - DefineProperBufferLength, 705
  - mElementOffsets, 706
  - mElementOffsets1, 706
  - mWriter, 707
  - mXMax, 707
  - mXMin, 707
  - mYMax, 707
  - mYMin, 707
  - mZMax, 707
  - mZMin, 707
  - mSPFile, 707
  - SetFile, 705
  - SetFileName, 705
  - SetStream, 705
  - StreamImageWriter, 704
  - Write, 705
  - WriteImageInformation, 706
  - WriteImageSubregionRAW, 706
  - WriteRawHeader, 706
- gdcmm::String
  - const\_iterator, 709
  - const\_reference, 709
  - const\_reverse\_iterator, 709
  - difference\_type, 709
  - IsValid, 710
  - iterator, 709
  - operator const char \*, 710
  - pointer, 709
  - reference, 709
  - reverse\_iterator, 709
  - size\_type, 709
  - String, 709, 710
  - Trim, 710
  - Truncate, 710
  - value\_type, 709
- gdcmm::String< TDelimiter, TMaxLength, TPadChar >, 707
- gdcmm::StringFilter, 710
  - ~StringFilter, 711
  - ExecuteQuery, 711
  - FromString, 712
  - GetFile, 712
  - SetDicts, 712
  - SetFile, 712
  - StringFilter, 711
  - ToString, 712
  - ToStringPair, 712
  - UseDictAlways, 713
- gdcmm::Study, 713
  - Study, 713
- gdcmm::Subject, 713
  - ~Subject, 715
  - AddObserver, 715
  - GetCommand, 715
  - HasObserver, 715
  - InvokeEvent, 715



- RemoveAllObservers, [715](#)
- RemoveObserver, [716](#)
- Subject, [715](#)
- gdcmm::Surface, [716](#)
  - ~Surface, [719](#)
  - GetAlgorithmFamily, [719](#)
  - GetAlgorithmName, [719](#)
  - GetAlgorithmVersion, [719](#)
  - GetAxisOfRotation, [719](#)
  - GetCenterOfRotation, [719](#)
  - GetFiniteVolume, [720](#)
  - GetManifold, [720](#)
  - GetMaximumPointDistance, [720](#)
  - GetMeanPointDistance, [720](#)
  - GetMeshPrimitive, [720](#)
  - GetNumberOfSurfacePoints, [720](#)
  - GetNumberOfVectors, [720](#)
  - GetPointCoordinatesData, [720](#)
  - GetPointPositionAccuracy, [720](#)
  - GetPointsBoundingBoxCoordinates, [720](#)
  - GetProcessingAlgorithm, [720](#)
  - GetRecommendedDisplayCIELabValue, [720](#)
  - GetRecommendedDisplayGrayscaleValue, [720](#)
  - GetRecommendedPresentationOpacity, [720](#)
  - GetRecommendedPresentationType, [720](#)
  - GetSTATES, [721](#)
  - GetSTATESString, [721](#)
  - GetSurfaceComments, [721](#)
  - GetSurfaceNumber, [721](#)
  - GetSurfaceProcessing, [721](#)
  - GetSurfaceProcessingDescription, [721](#)
  - GetSurfaceProcessingRatio, [721](#)
  - GetVIEWType, [721](#)
  - GetVIEWTypeString, [721](#)
  - GetVectorAccuracy, [721](#)
  - GetVectorCoordinateData, [721](#)
  - GetVectorDimensionality, [721](#)
  - NO, [719](#)
  - POINTS, [719](#)
  - STATES, [719](#)
  - STATES\_END, [719](#)
  - SURFACE, [719](#)
  - SetAlgorithmFamily, [721](#)
  - SetAlgorithmName, [721](#)
  - SetAlgorithmVersion, [721](#)
  - SetAxisOfRotation, [721](#)
  - SetCenterOfRotation, [721](#)
  - SetFiniteVolume, [721](#)
  - SetManifold, [721](#)
  - SetMaximumPointDistance, [721](#)
  - SetMeanPointDistance, [721](#)
  - SetMeshPrimitive, [721](#)
  - SetNumberOfSurfacePoints, [721](#)
  - SetNumberOfVectors, [721](#)
  - SetPointCoordinatesData, [721](#)
  - SetPointPositionAccuracy, [722](#)
  - SetPointsBoundingBoxCoordinates, [722](#)
  - SetProcessingAlgorithm, [722](#)
  - SetRecommendedDisplayCIELabValue, [722](#)
  - SetRecommendedDisplayGrayscaleValue, [722](#)
  - SetRecommendedPresentationOpacity, [722](#)
  - SetRecommendedPresentationType, [722](#)
  - SetSurfaceComments, [722](#)
  - SetSurfaceNumber, [722](#)
  - SetSurfaceProcessing, [722](#)
  - SetSurfaceProcessingDescription, [722](#)
  - SetSurfaceProcessingRatio, [722](#)
  - SetVectorAccuracy, [722](#)
  - SetVectorCoordinateData, [722](#)
  - SetVectorDimensionality, [722](#)
  - Surface, [719](#)
  - UNKNOWN, [719](#)
  - VIEWType, [719](#)
  - VIEWType\_END, [719](#)
  - WIREFRAME, [719](#)
  - YES, [719](#)
- gdcmm::SurfaceHelper, [722](#)
  - ColorArray, [723](#)
  - RGBToRecommendedDisplayCIELab, [724](#)
  - RGBToRecommendedDisplayGrayscale, [724](#)
  - RecommendedDisplayCIELabToRGB, [723](#)
- gdcmm::SurfaceReader, [725](#)
  - ~SurfaceReader, [727](#)
  - GetNumberOfSurfaces, [727](#)
  - Read, [727](#)
  - ReadPointMacro, [727](#)
  - ReadSurface, [727](#)
  - ReadSurfaces, [727](#)
  - SurfaceReader, [727](#)
- gdcmm::SurfaceWriter, [727](#)
  - ~SurfaceWriter, [729](#)
  - ComputeNumberOfSurfaces, [729](#)
  - GetNumberOfSurfaces, [729](#)
  - NumberOfSurfaces, [729](#)
  - PrepareWrite, [729](#)
  - PrepareWritePointMacro, [729](#)
  - SetNumberOfSurfaces, [729](#)
  - SurfaceWriter, [729](#)
  - Write, [729](#)
- gdcmm::SwapCode, [729](#)
  - BadBigEndian, [730](#)
  - BadLittleEndian, [730](#)
  - BigEndian, [730](#)
  - GetIndex, [731](#)
  - GetSwapCodeString, [731](#)
  - LittleEndian, [730](#)
  - operator SwapCode::SwapCodeType, [731](#)
  - operator<<, [731](#)

- SwapCode, 731
- SwapCodeType, 730
- Unknown, 730
- gdcmm::SwapperDoOp, 731
  - Swap, 731
  - SwapArray, 731
- gdcmm::SwapperNoOp, 732
  - Swap, 732
  - SwapArray, 732
- gdcmm::System, 732
  - DeleteDirectory, 733
  - EncodeBytes, 733
  - FileExists, 733
  - FileIsDirectory, 734
  - FileIsSymlink, 734
  - FileSize, 734
  - FileTime, 734
  - FormatDateTime, 734
  - GetCWD, 735
  - GetCurrentDateTime, 734
  - GetCurrentModuleFileName, 734
  - GetCurrentProcessFileName, 735
  - GetCurrentResourcesDirectory, 735
  - GetHostName, 735
  - GetLastError, 735
  - GetLocaleCharSet, 735
  - GetPermissions, 735
  - GetTimeZoneOffsetFromUTC, 735
  - MakeDirectory, 735
  - ParseDateTime, 735, 736
  - RemoveFile, 736
  - SetPermissions, 736
  - StrCaseCmp, 736
  - StrNCaseCmp, 736
  - StrSep, 736
  - StrTokR, 736
- gdcmm::Table, 736
  - ~Table, 737
  - GetTableEntry, 737
  - InsertEntry, 737
  - MapTableEntry, 737
  - operator<<, 737
  - Table, 737
- gdcmm::TableEntry, 737
  - ~TableEntry, 738
  - TableEntry, 738
- gdcmm::TableReader, 738
  - ~TableReader, 739
  - CharacterDataHandler, 739
  - EndElement, 739
  - GetDefs, 739
  - GetFilename, 739
  - HandleIOD, 739
  - HandleIODEntry, 739
  - HandleMacro, 739
  - HandleMacroEntry, 739
  - HandleMacroEntryDescription, 739
  - HandleModule, 740
  - HandleModuleEntry, 740
  - HandleModuleEntryDescription, 740
  - HandleModuleInclude, 740
  - Read, 740
  - SetFilename, 740
  - StartElement, 740
  - TableReader, 739
- gdcmm::Tag, 741
  - bytes, 748
  - GetElement, 744
  - GetElementTag, 744
  - GetGroup, 744
  - GetLength, 744
  - GetPrivateCreator, 744
  - IsGroupLength, 744
  - IsGroupXX, 744
  - IsIllegal, 745
  - IsPrivate, 745
  - IsPrivateCreator, 745
  - IsPublic, 745
  - operator!=, 745
  - operator<, 745
  - operator<<, 748
  - operator<=, 745
  - operator>>, 748
  - operator=, 745
  - operator==, 745
  - operator[], 746
  - PrintAsContinuousString, 746
  - PrintAsContinuousUpperCaseString, 746
  - PrintAsPipeSeparatedString, 746
  - Read, 746
  - ReadFromCommaSeparatedString, 746
  - ReadFromContinuousString, 746
  - ReadFromPipeSeparatedString, 746
  - SetElement, 747
  - SetElementTag, 747
  - SetGroup, 747
  - SetPrivateCreator, 747
  - Tag, 743
  - tag, 748
  - tags, 748
  - Write, 747
- gdcmm::TagPath, 748
  - ~TagPath, 749
  - ConstructFromString, 749
  - ConstructFromTagList, 749
  - IsValid, 749
  - Print, 749
  - Push, 749

- TagPath, [749](#)
- gdcmm::Testing, [749](#)
  - ~Testing, [751](#)
  - ComputeFileMD5, [751](#)
  - ComputeMD5, [751](#)
  - GetDataExtraRoot, [751](#)
  - GetDataRoot, [751](#)
  - GetFileName, [751](#)
  - GetFileNames, [751](#)
  - GetLossyFlagFromFile, [752](#)
  - GetMD5DataImage, [752](#)
  - GetMD5DataImages, [752](#)
  - GetMD5FromBrokenFile, [752](#)
  - GetMD5FromFile, [752](#)
  - GetMediaStorageDataFile, [752](#)
  - GetMediaStorageDataFiles, [752](#)
  - GetMediaStorageFromFile, [752](#)
  - GetNumberOfFileNames, [752](#)
  - GetNumberOfMD5DataImages, [752](#)
  - GetNumberOfMediaStorageDataFiles, [752](#)
  - GetPixelSpacingDataRoot, [752](#)
  - GetSelectedPrivateGroupOffsetFromFile, [752](#)
  - GetSelectedTagsOffsetFromFile, [753](#)
  - GetSourceDirectory, [753](#)
  - GetStreamOffsetFromFile, [753](#)
  - GetTempDirectory, [753](#)
  - GetTempDirectoryW, [753](#)
  - GetTempFilename, [753](#)
  - GetTempFilenameW, [753](#)
  - MD5DataImagesType, [750](#)
  - MediaStorageDataFilesType, [750](#)
  - Print, [753](#)
  - Testing, [751](#)
- gdcmm::Trace, [754](#)
  - ~Trace, [755](#)
  - DebugOff, [755](#)
  - DebugOn, [755](#)
  - ErrorOff, [755](#)
  - ErrorOn, [755](#)
  - GetDebugFlag, [755](#)
  - GetDebugStream, [755](#)
  - GetErrorFlag, [755](#)
  - GetErrorStream, [755](#)
  - GetStream, [756](#)
  - GetWarningFlag, [756](#)
  - GetWarningStream, [756](#)
  - SetDebug, [756](#)
  - SetDebugStream, [756](#)
  - SetError, [756](#)
  - SetErrorStream, [756](#)
  - SetStream, [756](#)
  - SetStreamToFile, [756](#)
  - SetWarning, [756](#)
  - SetWarningStream, [756](#)
- Trace, [755](#)
- WarningOff, [757](#)
- WarningOn, [757](#)
- gdcmm::TransferSyntax, [767](#)
  - CT\_private\_ELE, [760](#)
  - CanStoreLossy, [760](#)
  - DeflatedExplicitVRLittleEndian, [759](#)
  - Explicit, [759](#)
  - ExplicitVRBigEndian, [759](#)
  - ExplicitVRLittleEndian, [759](#)
  - GetNegociatedType, [760](#)
  - GetString, [760](#)
  - GetSwapCode, [760](#)
  - GetTSSString, [760](#)
  - GetTSType, [760](#)
  - Implicit, [759](#)
  - ImplicitVRBigEndianACRNEMA, [760](#)
  - ImplicitVRBigEndianPrivateGE, [759](#)
  - ImplicitVRLittleEndian, [759](#)
  - IsEncapsulated, [760](#)
  - IsEncoded, [760](#)
  - IsExplicit, [761](#)
  - IsImplicit, [761](#)
  - IsLossless, [761](#)
  - IsLossy, [761](#)
  - IsValid, [761](#)
  - JPEG2000, [760](#)
  - JPEG2000Lossless, [760](#)
  - JPEG2000Part2, [760](#)
  - JPEG2000Part2Lossless, [760](#)
  - JPEGBaselineProcess1, [759](#)
  - JPEGExtendedProcess2\_4, [759](#)
  - JPEGExtendedProcess3\_5, [759](#)
  - JPEGFullProgressionProcess10\_12, [759](#)
  - JPEGLSLossless, [759](#)
  - JPEGLSNearLossless, [759](#)
  - JPEGLosslessProcess14, [759](#)
  - JPEGLosslessProcess14\_1, [759](#)
  - JPEGSpectralSelectionProcess6\_8, [759](#)
  - JPIPRReferenced, [760](#)
  - MPEG2MainProfile, [760](#)
  - NegociatedType, [759](#)
  - operator TSType, [761](#)
  - operator<<, [761](#)
  - RLELossless, [760](#)
  - TS\_END, [760](#)
  - TSType, [759](#)
  - TransferSyntax, [760](#)
  - Unknown, [759](#)
- gdcmm::Type, [764](#)
  - GetTypeString, [765](#)
  - GetTypeType, [765](#)
  - operator TypeType, [765](#)
  - operator<<, [765](#)

- T1, [765](#)
- T1C, [765](#)
- T2, [765](#)
- T2C, [765](#)
- T3, [765](#)
- Type, [765](#)
- TypeType, [765](#)
- UNKNOWN, [765](#)
- gdcmm::UI, [765](#)
  - Internal, [766](#)
  - operator<<, [766](#)
- gdcmm::UIDGenerator, [766](#)
  - Generate, [767](#)
  - GenerateUUID, [767](#)
  - GetGDCMUID, [767](#)
  - GetRoot, [767](#)
  - IsValid, [767](#)
  - SetRoot, [767](#)
  - UIDGenerator, [767](#)
- gdcmm::UIDs, [768](#)
  - AmbulatoryECGWaveformStorage, [775](#)
  - AudioSRStorageTrialRetired, [776](#)
  - BasicAnnotationBoxSOPClass, [774](#)
  - BasicColorImageBoxSOPClass, [774](#)
  - BasicColorPrintManagementMetaSOPClass, [774](#)
  - BasicFilmBoxSOPClass, [774](#)
  - BasicFilmSessionSOPClass, [774](#)
  - BasicGrayscaleImageBoxSOPClass, [774](#)
  - BasicGrayscalePrintManagementMetaSOPClass, [774](#)
  - BasicPrintImageOverlayBoxSOPClassRetired, [775](#)
  - BasicStudyContentNotificationSOPClassRetired, [774](#)
  - BasicTextSRStorage, [776](#)
  - BasicVoiceAudioWaveformStorage, [775](#)
  - BlendingSoftcopyPresentationStateStorageSOPClass, [776](#)
  - BreastImagingRelevantPatientInformationQuery, [777](#)
  - BreastTomosynthesisImageStorage, [779](#)
  - CTImageStorage, [775](#)
  - CardiacElectrophysiologyWaveformStorage, [775](#)
  - CardiacRelevantPatientInformationQuery, [778](#)
  - ChestCADSRStorage, [777](#)
  - ColorSoftcopyPresentationStateStorageSOPClass, [776](#)
  - ComprehensiveSRStorage, [776](#)
  - ComprehensiveSRStorageTrialRetired, [776](#)
  - ComputedRadiographyImageStorage, [775](#)
  - DICOMApplicationContextName, [774](#)
  - DICOMControlledTerminology, [774](#)
  - DICOMUIDRegistry, [774](#)
  - DeflatedExplicitVRLittleEndian, [772](#)
  - DeformableSpatialRegistrationStorage, [776](#)
  - DetachedInterpretationManagementSOPClassRetired, [774](#)
  - DetachedPatientManagementMetaSOPClassRetired, [774](#)
  - DetachedPatientManagementSOPClassRetired, [774](#)
  - DetachedResultsManagementMetaSOPClassRetired, [774](#)
  - DetachedResultsManagementSOPClassRetired, [774](#)
  - DetachedStudyManagementMetaSOPClassRetired, [774](#)
  - DetachedStudyManagementSOPClassRetired, [774](#)
  - DetachedVisitManagementSOPClassRetired, [774](#)
  - DetailSRStorageTrialRetired, [776](#)
  - dicomAETitle, [778](#)
  - dicomApplicationCluster, [778](#)
  - dicomAssociationAcceptor, [778](#)
  - dicomAssociationInitiator, [778](#)
  - dicomAuthorizedNodeCertificateReference, [778](#)
  - dicomConfigurationRoot, [778](#)
  - dicomDescription, [778](#)
  - dicomDevice, [778](#)
  - dicomDeviceName, [778](#)
  - dicomDeviceSerialNumber, [778](#)
  - dicomDevicesRoot, [778](#)
  - dicomHostname, [778](#)
  - dicomInstalled, [778](#)
  - dicomInstitutionAddress, [778](#)
  - dicomInstitutionDepartmentName, [778](#)
  - dicomInstitutionName, [778](#)
  - dicomIssuerOfPatientID, [778](#)
  - dicomManufacturer, [778](#)
  - dicomManufacturerModelName, [778](#)
  - dicomNetworkAE, [778](#)
  - dicomNetworkConnection, [779](#)
  - dicomNetworkConnectionReference, [778](#)
  - dicomPort, [778](#)
  - dicomPreferredCalledAETitle, [778](#)
  - dicomPreferredCallingAETitle, [778](#)
  - dicomPrimaryDeviceType, [778](#)
  - dicomRelatedDeviceReference, [778](#)
  - dicomSOPClass, [778](#)
  - dicomSoftwareVersion, [778](#)
  - dicomStationName, [778](#)
  - dicomSupportedCharacterSet, [778](#)
  - dicomTLSCyphersuite, [778](#)
  - dicomThisNodeCertificateReference, [778](#)
  - dicomTransferCapability, [779](#)
  - dicomTransferRole, [778](#)
  - dicomTransferSyntax, [778](#)
  - dicomUniqueAETitle, [779](#)
  - dicomUniqueAETitlesRegistryRoot, [778](#)
  - dicomVendorData, [778](#)

- DigitalIntraoralXRayImageStorageForPresentation, [775](#)
- DigitalIntraoralXRayImageStorageForProcessing, [775](#)
- DigitalMammographyXRayImageStorageFor↔ Presentation, [775](#)
- DigitalMammographyXRayImageStorageFor↔ Processing, [775](#)
- DigitalXRayImageStorageForPresentation, [775](#)
- DigitalXRayImageStorageForProcessing, [775](#)
- EncapsulatedCDASStorage, [777](#)
- EncapsulatedPDFStorage, [777](#)
- EnhancedCTImageStorage, [775](#)
- EnhancedMRIImageStorage, [775](#)
- EnhancedSRStorage, [776](#)
- EnhancedUSVolumeStorage, [779](#)
- EnhancedXAImageStorage, [776](#)
- EnhancedXRFImageStorage, [776](#)
- ExplicitVRBigEndian, [772](#)
- ExplicitVRLittleEndian, [772](#)
- GeneralECGWaveformStorage, [775](#)
- GeneralPurposePerformedProcedureStepSOP↔ Class, [777](#)
- GeneralPurposeScheduledProcedureStepSOP↔ Class, [777](#)
- GeneralPurposeWorklistInformationModelFIND, [777](#)
- GeneralPurposeWorklistManagementMetaSOP↔ Class, [777](#)
- GeneralRelevantPatientInformationQuery, [777](#)
- GetName, [785](#)
- GetNumberOfTransferSyntaxStrings, [785](#)
- GetString, [786](#)
- GetTransferSyntaxString, [786](#)
- GetTransferSyntaxStrings, [786](#)
- GetUIDName, [786](#)
- GetUIDString, [786](#)
- GrayscaleSoftcopyPresentationStateStorageSOP↔ Class, [776](#)
- HangingProtocolInformationModelFIND, [778](#)
- HangingProtocolInformationModelMOVE, [778](#)
- HangingProtocolStorage, [778](#)
- HardcopyColorImageStorageSOPClassRetired, [775](#)
- HardcopyGrayscaleImageStorageSOPClassRetired, [775](#)
- HemodynamicWaveformStorage, [775](#)
- ICBM452T1FrameofReference, [774](#)
- ICBMSingleSubjectMRIFrameofReference, [774](#)
- ImageOverlayBoxSOPClassRetired, [775](#)
- ImplicitVRLittleEndianDefaultTransferSyntaxforDIC↔ OM, [772](#)
- InstanceAvailabilityNotificationSOPClass, [777](#)
- JPEG2000ImageCompression, [773](#)
- JPEG2000ImageCompressionLosslessOnly, [773](#)
- JPEG2000Part2MulticomponentImageCompression, [773](#)
- JPEG2000Part2MulticomponentImageCompression↔ LosslessOnly, [773](#)
- JPEGBaselineProcess1DefaultTransferSyntaxfor↔ LossyJPEG8BitImageCompression, [772](#)
- JPEGExtendedHierarchicalProcess1618Retired, [773](#)
- JPEGExtendedHierarchicalProcess1719Retired, [773](#)
- JPEGExtendedProcess24DefaultTransferSyntaxfor↔ LossyJPEG12BitImageCompressionProcess4only, [772](#)
- JPEGExtendedProcess35Retired, [772](#)
- JPEGFullProgressionHierarchicalProcess2426↔ Retired, [773](#)
- JPEGFullProgressionHierarchicalProcess2527↔ Retired, [773](#)
- JPEGFullProgressionNonHierarchicalProcess1012↔ Retired, [772](#)
- JPEGFullProgressionNonHierarchicalProcess1113↔ Retired, [772](#)
- JPEGLSLosslessImageCompression, [773](#)
- JPEGLSLossyNearLosslessImageCompression, [773](#)
- JPEGLosslessHierarchicalProcess28Retired, [773](#)
- JPEGLosslessHierarchicalProcess29Retired, [773](#)
- JPEGLosslessNonHierarchicalFirstOrderPrediction↔ Process14SelectionValue1DefaultTransfer↔ SyntaxforLosslessJPEGImageCompression, [773](#)
- JPEGLosslessNonHierarchicalProcess14, [772](#)
- JPEGLosslessNonHierarchicalProcess15Retired, [773](#)
- JPEGSpectralSelectionHierarchicalProcess2022↔ Retired, [773](#)
- JPEGSpectralSelectionHierarchicalProcess2123↔ Retired, [773](#)
- JPEGSpectralSelectionNonHierarchicalProcess68↔ Retired, [772](#)
- JPEGSpectralSelectionNonHierarchicalProcess79↔ Retired, [772](#)
- JPIPRReferenced, [773](#)
- JPIPRReferencedDeflate, [773](#)
- KeyObjectSelectionDocumentStorage, [777](#)
- MPEG2MainProfileMainLevel, [773](#)
- MRIImageStorage, [775](#)
- MRSpectroscopyStorage, [775](#)
- MammographyCADSRStorage, [776](#)
- MediaCreationManagementSOPClassUID, [775](#)
- MediaStorageDirectoryStorage, [773](#)
- ModalityPerformedProcedureStepNotificationSOP↔ Class, [774](#)
- ModalityPerformedProcedureStepRetrieveSOP↔ Class, [774](#)
- ModalityPerformedProcedureStepSOPClass, [774](#)
- ModalityWorklistInformationModelFIND, [777](#)
- MultiframeGrayscaleByteSecondaryCaptureImage↔

- Storage, [775](#)
- MultiframeGrayscaleWordSecondaryCapture↔  
ImageStorage, [775](#)
- MultiframeSingleBitSecondaryCaptureImage↔  
Storage, [775](#)
- MultiframeTrueColorSecondaryCaptureImage↔  
Storage, [775](#)
- NuclearMedicineImageStorage, [776](#)
- NuclearMedicineImageStorageRetired, [775](#)
- operator TSType, [786](#)
- OphthalmicPhotography16BitImageStorage, [776](#)
- OphthalmicPhotography8BitImageStorage, [776](#)
- OphthalmicTomographyImageStorage, [776](#)
- PatientRootQueryRetrieveInformationModelFIND,  
[777](#)
- PatientRootQueryRetrieveInformationModelGE↔  
T, [777](#)
- PatientRootQueryRetrieveInformationModelMOVE,  
[777](#)
- PatientStudyOnlyQueryRetrieveInformationModel↔  
FINDRetired, [777](#)
- PatientStudyOnlyQueryRetrieveInformationModel↔  
GETRetired, [777](#)
- PatientStudyOnlyQueryRetrieveInformationModel↔  
MOVERetired, [777](#)
- PositronEmissionTomographyImageStorage, [777](#)
- PresentationLUTSOPClass, [775](#)
- PrintJobSOPClass, [774](#)
- PrintQueueManagementSOPClassRetired, [775](#)
- PrintQueueSOPInstanceRetired, [775](#)
- PrinterConfigurationRetrieveSOPClass, [774](#)
- PrinterConfigurationRetrieveSOPInstance, [774](#)
- PrinterSOPClass, [774](#)
- PrinterSOPInstance, [774](#)
- ProceduralEventLoggingSOPClass, [774](#)
- ProceduralEventLoggingSOPInstance, [774](#)
- ProcedureLogStorage, [776](#)
- ProductCharacteristicsQuerySOPClass, [778](#)
- PseudoColorSoftcopyPresentationStateStorageS↔  
OPClass, [776](#)
- PullPrintRequestSOPClassRetired, [775](#)
- PullStoredPrintManagementMetaSOPClassRetired,  
[775](#)
- RFC2557MIMEencapsulation, [773](#)
- RLELossless, [773](#)
- RTBeamsDeliveryInstructionStorageSupplement74↔  
FrozenDraft, [777](#)
- RTBeamsTreatmentRecordStorage, [777](#)
- RTBrachyTreatmentRecordStorage, [777](#)
- RTConventionalMachineVerificationSupplement74↔  
FrozenDraft, [777](#)
- RTDoseStorage, [777](#)
- RTImageStorage, [777](#)
- RTIonBeamsTreatmentRecordStorage, [777](#)
- RTIonMachineVerificationSupplement74FrozenDraft,  
[777](#)
- RTIonPlanStorage, [777](#)
- RTPlanStorage, [777](#)
- RTStructureSetStorage, [777](#)
- RTTreatmentSummaryRecordStorage, [777](#)
- RawDataStorage, [776](#)
- RealWorldValueMappingStorage, [776](#)
- ReferencedColorPrintManagementMetaSOPClass↔  
Retired, [774](#)
- ReferencedGrayscalePrintManagementMetaSOP↔  
ClassRetired, [774](#)
- ReferencedImageBoxSOPClassRetired, [774](#)
- SPM2AVG152PDFFrameofReference, [773](#)
- SPM2AVG152T1FrameofReference, [773](#)
- SPM2AVG152T2FrameofReference, [773](#)
- SPM2AVG305T1FrameofReference, [773](#)
- SPM2BRAINMASKFrameofReference, [773](#)
- SPM2CSFFFrameofReference, [773](#)
- SPM2EPIFrameofReference, [773](#)
- SPM2FILT1FrameofReference, [773](#)
- SPM2GRAYFrameofReference, [773](#)
- SPM2PDFFrameofReference, [773](#)
- SPM2PETFrameofReference, [773](#)
- SPM2SINGLESUBJT1FrameofReference, [773](#)
- SPM2SPECTFrameofReference, [773](#)
- SPM2T1FrameofReference, [773](#)
- SPM2T2FrameofReference, [773](#)
- SPM2TRANSMFrameofReference, [773](#)
- SPM2WHITEFrameofReference, [773](#)
- SecondaryCaptureImageStorage, [775](#)
- SegmentationStorage, [776](#)
- SetFromUID, [786](#)
- SpatialFiducialsStorage, [776](#)
- SpatialRegistrationStorage, [776](#)
- StandaloneCurveStorageRetired, [775](#)
- StandaloneModalityLUTStorageRetired, [776](#)
- StandaloneOverlayStorageRetired, [775](#)
- StandalonePETCurveStorageRetired, [777](#)
- StandaloneVOILUTStorageRetired, [776](#)
- StereometricRelationshipStorage, [776](#)
- StorageCommitmentPullModelSOPClassRetired,  
[774](#)
- StorageCommitmentPullModelSOPInstanceRetired,  
[774](#)
- StorageCommitmentPushModelSOPClass, [774](#)
- StorageCommitmentPushModelSOPInstance, [774](#)
- StorageServiceClass, [774](#)
- StoredPrintStorageSOPClassRetired, [775](#)
- StudyComponentManagementSOPClassRetired,  
[774](#)
- StudyRootQueryRetrieveInformationModelIFIND, [777](#)
- StudyRootQueryRetrieveInformationModelIGET, [777](#)



StudyRootQueryRetrieveInformationModelMOVE, [777](#)  
SubstanceAdministrationLoggingSOPClass, [774](#)  
SubstanceAdministrationLoggingSOPInstance, [774](#)  
SubstanceApprovalQuerySOPClass, [778](#)  
SurfaceSegmentationStorage, [779](#)  
TSName, [772](#)  
TSType, [779](#)  
TalairachBrainAtlasFrameofReference, [773](#)  
TextSRStorageTrialRetired, [776](#)  
TransferSyntaxStringsType, [772](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_1, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_10, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_11, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_12, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_13, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_14, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_15, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_16, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_17, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_18, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_19, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_2, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_20, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_21, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_22, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_23, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_24, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_25, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_26, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_27, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_28, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_29, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_3, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_30, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_31, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_4, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_5, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_6, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_7, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_8, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_9, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_1, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_2, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_3, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_4, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_5, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_6, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_7, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_8, [785](#)  
uid\_1\_2\_840\_10008\_1\_1, [779](#)  
uid\_1\_2\_840\_10008\_1\_2, [779](#)  
uid\_1\_2\_840\_10008\_1\_20\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_1\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_2, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_2\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_1, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_1\_99, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_2, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_100, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_50, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_51, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_52, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_53, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_54, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_55, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_56, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_57, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_58, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_59, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_60, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_61, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_62, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_63, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_64, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_65, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_66, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_70, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_80, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_81, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_90, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_91, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_92, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_93, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_94, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_95, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_5, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_6\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_6\_2, [780](#)  
uid\_1\_2\_840\_10008\_1\_3\_10, [780](#)  
uid\_1\_2\_840\_10008\_1\_40, [780](#)  
uid\_1\_2\_840\_10008\_1\_40\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_42, [780](#)  
uid\_1\_2\_840\_10008\_1\_42\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_1, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_10, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_11, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_12, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_13, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_14, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_15, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_16, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_17, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_18, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_2, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_3, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_4, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_5, [780](#)

uid\_1\_2\_840\_10008\_1\_4\_1\_6, 780  
 uid\_1\_2\_840\_10008\_1\_4\_1\_7, 780  
 uid\_1\_2\_840\_10008\_1\_4\_1\_8, 780  
 uid\_1\_2\_840\_10008\_1\_4\_1\_9, 780  
 uid\_1\_2\_840\_10008\_1\_4\_2\_1, 780  
 uid\_1\_2\_840\_10008\_1\_4\_2\_2, 780  
 uid\_1\_2\_840\_10008\_1\_9, 780  
 uid\_1\_2\_840\_10008\_2\_16\_4, 780  
 uid\_1\_2\_840\_10008\_2\_6\_1, 780  
 uid\_1\_2\_840\_10008\_3\_1\_1\_1, 780  
 uid\_1\_2\_840\_10008\_3\_1\_2\_1\_1, 780  
 uid\_1\_2\_840\_10008\_3\_1\_2\_1\_4, 780  
 uid\_1\_2\_840\_10008\_3\_1\_2\_2\_1, 780  
 uid\_1\_2\_840\_10008\_3\_1\_2\_3\_1, 780  
 uid\_1\_2\_840\_10008\_3\_1\_2\_3\_2, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_3\_3, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_3\_4, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_3\_5, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_5\_1, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_5\_4, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_5\_5, 781  
 uid\_1\_2\_840\_10008\_3\_1\_2\_6\_1, 781  
 uid\_1\_2\_840\_10008\_4\_2, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_14, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_15, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_16, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_16\_376, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_17, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_17\_376, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_18, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_18\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_2, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_22, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_23, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_24, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_24\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_25, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_26, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_27, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_29, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_30, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_31, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_32, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_33, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_4, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_4\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_4\_2, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_9, 781  
 uid\_1\_2\_840\_10008\_5\_1\_1\_9\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_10, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_1, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_2, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_3, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_4, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_128, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_129, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_3, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_3, 785  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2\_1, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3, 781  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_20, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_1, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_2, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_3, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_4, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_5, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_6, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_7, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_8, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_9, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_5, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_2, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_3, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_4, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_5, 785  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_67, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6\_1, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6\_2, 785  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7, 782  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_1, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_1\_1, 783  
 uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2, 783



- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2\_1, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_3, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4\_1, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_1, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_2, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_3, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_4, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_6, [785](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_2, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_2, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_3, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_4, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_8, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_1, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_11, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_2, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_22, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_3, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_33, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_4, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_40, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_50, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_59, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_65, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_67, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_2, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_3, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_2\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_3\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_4\_1, [782](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_1, [783](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_31, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_32, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_32\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_32\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_32\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_33, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_4, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_34\_5, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_37\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_37\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_37\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_38\_1, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_38\_2, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_38\_3, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_41, [784](#)
- uid\_1\_2\_840\_10008\_5\_1\_4\_42, [784](#)
- UltrasoundImageStorage, [775](#)
- UltrasoundImageStorageRetired, [775](#)
- UltrasoundMultiframeImageStorage, [775](#)
- UltrasoundMultiframeImageStorageRetired, [775](#)
- UnifiedProcedureStepEventSOPClass, [777](#)
- UnifiedProcedureStepPullSOPClass, [777](#)
- UnifiedProcedureStepPushSOPClass, [777](#)
- UnifiedProcedureStepWatchSOPClass, [777](#)
- UnifiedWorklistandProcedureStepSOPInstance, [777](#)
- UnifiedWorklistandProcedureStepServiceClass, [777](#)
- VLEndoscopicImageStorage, [776](#)
- VImageStorageTrialRetired, [776](#)
- VLMicroscopicImageStorage, [776](#)
- VLMultiframeImageStorageTrialRetired, [776](#)
- VLPhotographicImageStorage, [776](#)
- VLSlideCoordinatesMicroscopicImageStorage, [776](#)
- VLWholeSlideMicroscopyImageStorage, [779](#)
- VOILUTBoxSOPClass, [775](#)
- VerificationSOPClass, [772](#)
- VideoEndoscopicImageStorage, [776](#)
- VideoMicroscopicImageStorage, [776](#)
- VideoPhotographicImageStorage, [776](#)
- WaveformStorageTrialRetired, [775](#)
- XMLEncoding, [773](#)
- XRay3DAngiographicImageStorage, [776](#)
- XRay3DCraniofacialImageStorage, [776](#)
- XRayAngiographicBiPlaneImageStorageRetired, [776](#)
- XRayAngiographicImageStorage, [776](#)
- XRayRadiationDoseSRStorage, [777](#)
- XRayRadiofluoroscopicImageStorage, [776](#)
- gdcmm::UNExplicitDataElement, [830](#)
  - GetLength, [831](#)
  - Read, [831](#)
  - ReadPreValue, [831](#)
  - ReadValue, [831](#)
  - ReadWithLength, [831](#)
- gdcmm::UNExplicitImplicitDataElement, [831](#)
  - GetLength, [833](#)
  - Read, [833](#)
  - ReadPreValue, [833](#)
  - ReadValue, [833](#)

- gdcmm::UUIDGenerator, 838
  - Generate, 838
  - IsValid, 838
- gdcmm::Unpacker12Bits, 833
  - Pack, 834
  - Unpack, 834
- gdcmm::Usage, 834
  - Conditional, 835
  - GetUsageString, 835
  - GetUsageType, 835
  - Invalid, 835
  - Mandatory, 835
  - operator UsageType, 835
  - operator<<, 835
  - Usage, 835
  - UsageType, 835
  - UserOption, 835
- gdcmm::UserEvent, 836
- gdcmm::VL, 844
  - GetLength, 845
  - GetVL16Max, 845
  - GetVL32Max, 845
  - IsOdd, 845
  - IsUndefined, 845
  - operator uint32\_t, 846
  - operator<<, 846
  - operator++, 846
  - operator+&, 846
  - Read, 846
  - Read16, 846
  - SetToUndefined, 846
  - Type, 845
  - VL, 845
  - Write, 846
  - Write16, 846
- gdcmm::VM, 846
  - Compatible, 849
  - GetIndex, 849
  - GetLength, 849
  - GetNumberOfElementsFromArray, 849
  - GetVMString, 849
  - GetVMType, 849
  - GetVMTypeFromLength, 850
  - IsValid, 850
  - operator VMType, 850
  - operator<<, 850
  - VM, 849
  - VM0, 848
  - VM1, 848
  - VM10, 848
  - VM12, 848
  - VM16, 848
  - VM18, 848
  - VM1\_2, 849
  - VM1\_3, 849
  - VM1\_32, 849
  - VM1\_4, 849
  - VM1\_5, 849
  - VM1\_8, 849
  - VM1\_99, 849
  - VM1\_n, 849
  - VM2, 848
  - VM24, 848
  - VM256, 849
  - VM28, 848
  - VM2\_2n, 849
  - VM2\_n, 849
  - VM3, 848
  - VM30\_30n, 849
  - VM32, 848
  - VM35, 848
  - VM3\_3n, 849
  - VM3\_4, 849
  - VM3\_n, 849
  - VM4, 848
  - VM47\_47n, 849
  - VM4\_4n, 849
  - VM5, 848
  - VM6, 848
  - VM6\_6n, 849
  - VM7\_7n, 849
  - VM8, 848
  - VM9, 848
  - VM99, 849
  - VM\_END, 849
  - VMType, 848
- gdcmm::VMToLength< T >, 850
- gdcmm::VR, 850
  - AE, 852
  - AS, 852
  - AT, 852
  - CS, 852
  - CanDisplay, 854
  - Compatible, 854
  - DA, 852
  - DS, 853
  - DT, 853
  - FD, 853
  - FL, 853
  - GetLength, 854
  - GetSize, 854
  - GetSizeof, 854
  - GetVRString, 854
  - GetVRStringFromFile, 854
  - GetVRType, 854
  - GetVRTypeFromFile, 854
  - INVALID, 852
  - IS, 853

- IsASCII, [854](#)
- IsASCII2, [854](#)
- IsBinary, [854](#)
- IsBinary2, [854](#)
- IsDual, [854](#)
- IsSwap, [854](#)
- IsVRFile, [854](#)
- IsValid, [854](#)
- LO, [853](#)
- LT, [853](#)
- OB, [853](#)
- OB\_OW, [853](#)
- OD, [853](#)
- OF, [853](#)
- OW, [853](#)
- operator VRType, [854](#)
- operator<=, [855](#)
- PN, [853](#)
- Read, [855](#)
- SH, [853](#)
- SL, [853](#)
- SQ, [853](#)
- SS, [853](#)
- ST, [853](#)
- TM, [853](#)
- UI, [853](#)
- UL, [853](#)
- UN, [853](#)
- US, [853](#)
- US\_SS, [853](#)
- US\_SS\_OW, [853](#)
- UT, [853](#)
- VL16, [853](#)
- VL32, [853](#)
- VR, [854](#)
- VR\_END, [853](#)
- VR\_VM1, [853](#)
- VRALL, [853](#)
- VRASCII, [853](#)
- VRBINARY, [853](#)
- VRType, [852](#)
- Write, [855](#)
- gdcmm::VR16ExplicitDataElement, [855](#)
  - GetLength, [856](#)
  - Read, [856](#)
  - ReadPreValue, [856](#)
  - ReadValue, [857](#)
  - ReadWithLength, [857](#)
- gdcmm::VRToEncoding< T >, [857](#)
- gdcmm::VRToType< T >, [857](#)
- gdcmm::VRVLSIZE< 0 >, [858](#)
  - Read, [858](#)
  - Write, [858](#)
- gdcmm::VRVLSIZE< 1 >, [858](#)
  - Read, [858](#)
  - Write, [858](#)
- gdcmm::VRVLSIZE< T >, [858](#)
- gdcmm::Validate, [839](#)
  - ~Validate, [839](#)
  - F, [840](#)
  - GetValidatedFile, [840](#)
  - SetFile, [840](#)
  - V, [840](#)
  - Validate, [839](#)
  - Validation, [840](#)
- gdcmm::Value, [840](#)
  - ~Value, [841](#)
  - Clear, [841](#)
  - DataElement, [842](#)
  - GetLength, [842](#)
  - operator==, [842](#)
  - SetLength, [842](#)
  - SetLengthOnly, [842](#)
  - Value, [841](#)
- gdcmm::ValueIO
  - Read, [843](#)
  - Write, [843](#)
- gdcmm::ValueIO< TDE, TSwap, TType >, [842](#)
- gdcmm::Version, [843](#)
  - ~Version, [843](#)
  - GetBuildVersion, [843](#)
  - GetMajorVersion, [844](#)
  - GetMinorVersion, [844](#)
  - GetVersion, [844](#)
  - operator<=, [844](#)
  - Print, [844](#)
  - Version, [843](#)
- gdcmm::Waveform, [916](#)
  - Waveform, [916](#)
- gdcmm::Writer, [916](#)
  - ~Writer, [920](#)
  - CheckFileMetaInformationOff, [920](#)
  - CheckFileMetaInformationOn, [920](#)
  - GetFile, [920](#)
  - GetStreamPtr, [920](#)
  - Ofstream, [921](#)
  - SetCheckFileMetaInformation, [920](#)
  - SetFile, [920](#)
  - SetFileName, [920](#)
  - SetStream, [921](#)
  - SetWriteDataSetOnly, [921](#)
  - Stream, [921](#)
  - StreamImageWriter, [921](#)
  - Write, [921](#)
  - Writer, [920](#)
- gdcmm::XMLDictReader, [921](#)
  - ~XMLDictReader, [923](#)
  - CharacterDataHandler, [923](#)

- EndElement, [923](#)
- GetDict, [923](#)
- HandleDescription, [923](#)
- HandleEntry, [923](#)
- StartElement, [923](#)
- XMLDictReader, [923](#)
- gdcmm::XMLPrinter, [923](#)
  - ~XMLPrinter, [925](#)
  - F, [925](#)
  - GetPrintStyle, [925](#)
  - HandleBulkData, [925](#)
  - LOADBULKDATA, [925](#)
  - OnlyUUID, [925](#)
  - Print, [925](#)
  - PrintDataElement, [925](#)
  - PrintDataSet, [925](#)
  - PrintSQ, [925](#)
  - PrintStyle, [925](#)
  - PrintStyles, [925](#)
  - SetFile, [925](#)
  - SetStyle, [925](#)
  - XMLPrinter, [925](#)
- gdcmm::XMLPrivateDictReader, [926](#)
  - ~XMLPrivateDictReader, [927](#)
  - CharacterDataHandler, [927](#)
  - EndElement, [927](#)
  - GetPrivateDict, [927](#)
  - HandleDescription, [927](#)
  - HandleEntry, [927](#)
  - StartElement, [927](#)
  - XMLPrivateDictReader, [927](#)
- gdcmm::ignore\_char, [412](#)
  - ignore\_char, [413](#)
  - m\_char, [413](#)
- gdcmm::network, [131](#)
  - cMaxEventID, [136](#)
  - cMaxStateID, [136](#)
  - eAABORTPDUReceivedOpen, [135](#)
  - eAABORTRequest, [135](#)
  - eAASSOCIATE\_RQPDUreceived, [135](#)
  - eAASSOCIATERequestLocalUser, [135](#)
  - eAASSOCIATEResponseAccept, [135](#)
  - eAASSOCIATEResponseReject, [135](#)
  - eARELEASE\_RPPDUReceived, [135](#)
  - eARELEASE\_RQPDUReceivedOpen, [135](#)
  - eARELEASERequest, [135](#)
  - eARELEASEResponse, [135](#)
  - eARTIMTimerExpired, [135](#)
  - eASSOCIATE\_ACPDUreceived, [135](#)
  - eASSOCIATE\_RJPDUreceived, [135](#)
  - eEventDoesNotExist, [135](#)
  - EEventID, [135](#)
  - ePDATATFPDU, [135](#)
  - ePDATArequest, [135](#)
  - eSta10ReleaseCollisionAc, [136](#)
  - eSta11ReleaseCollisionRq, [136](#)
  - eSta12ReleaseCollisionAcLocal, [136](#)
  - eSta13AwaitingClose, [136](#)
  - eSta1Idle, [135](#)
  - eSta2Open, [135](#)
  - eSta3WaitLocalAssoc, [135](#)
  - eSta4LocalAssocDone, [135](#)
  - eSta5WaitRemoteAssoc, [135](#)
  - eSta6TransferReady, [135](#)
  - eSta7WaitRelease, [135](#)
  - eSta8WaitLocalRelease, [136](#)
  - eSta9ReleaseCollisionRqLocal, [136](#)
  - eStaDoesNotExist, [135](#)
  - EStateID, [135](#)
  - eTransportConnConfirmLocal, [135](#)
  - eTransportConnIndicLocal, [135](#)
  - eTransportConnectionClosed, [135](#)
  - eUnrecognizedPDUReceived, [135](#)
  - GetStateIndex, [136](#)
- gdcmm::network::AAAbortPDU, [139](#)
  - AAAbortPDU, [140](#)
  - IsLastFragment, [140](#)
  - Print, [140](#)
  - Read, [140](#)
  - SetReason, [140](#)
  - SetSource, [141](#)
  - Size, [141](#)
  - Write, [141](#)
- gdcmm::network::AAssociateACPDU, [141](#)
  - AAssociateACPDU, [143](#)
  - AAssociateRQPDU, [143](#)
  - AddPresentationContextAC, [143](#)
  - GetNumberOfPresentationContextAC, [143](#)
  - GetPresentationContextAC, [143](#)
  - GetUserInfo, [143](#)
  - InitFromRQ, [143](#)
  - IsLastFragment, [143](#)
  - Print, [143](#)
  - Read, [143](#)
  - SetCalledAETitle, [143](#)
  - SetCallingAETitle, [143](#)
  - Size, [143](#)
  - SizeType, [143](#)
  - Write, [143](#)
- gdcmm::network::AAssociateRJPDU, [144](#)
  - AAssociateRJPDU, [145](#)
  - IsLastFragment, [145](#)
  - Print, [145](#)
  - Read, [145](#)
  - Size, [145](#)
  - Write, [145](#)
- gdcmm::network::AAssociateRQPDU, [145](#)
  - AAssociateACPDU, [149](#)

- AAAssociateRQPDU, 147
- AddPresentationContext, 147
- GetCalledAETitle, 147
- GetCallingAETitle, 148
- GetNumberOfPresentationContext, 148
- GetPresentationContext, 148
- GetPresentationContextByAbstractSyntax, 148
- GetPresentationContextByID, 148
- GetPresentationContexts, 148
- GetReserved43\_74, 148
- GetUserInformation, 148
- IsAETitleValid, 148
- IsLastFragment, 148
- PresentationContextArrayType, 147
- Print, 148
- Read, 148
- SetCalledAETitle, 148
- SetCallingAETitle, 148
- SetUserInformation, 148
- Size, 149
- SizeType, 147
- Write, 149
- gdcmm::network::ARTIMTimer, 166
  - ARTIMTimer, 167
  - GetElapsedTime, 167
  - GetHasExpired, 167
  - GetTimeout, 167
  - SetTimeout, 167
  - Start, 167
  - Stop, 167
- gdcmm::network::AReleaseRPPDU, 163
  - AReleaseRPPDU, 164
  - IsLastFragment, 164
  - Print, 164
  - Read, 164
  - Size, 164
  - Write, 164
- gdcmm::network::AReleaseRQPDU, 164
  - AReleaseRQPDU, 165
  - IsLastFragment, 166
  - Print, 166
  - Read, 166
  - Size, 166
  - Write, 166
- gdcmm::network::AbstractSyntax, 150
  - AbstractSyntax, 151
  - GetAsDataElement, 151
  - GetName, 151
  - operator==, 151
  - Print, 151
  - Read, 151
  - SetName, 151
  - SetNameFromUID, 151
  - Size, 151
  - Write, 151
- gdcmm::network::ApplicationContext, 160
  - ApplicationContext, 160
  - GetName, 161
  - Print, 161
  - Read, 161
  - SetName, 161
  - Size, 161
  - Write, 161
- gdcmm::network::AsynchronousOperationsWindowSub, 168
  - AsynchronousOperationsWindowSub, 168
  - Print, 168
  - Read, 169
  - Size, 169
  - Write, 169
- gdcmm::network::BaseCompositeMessage, 196
  - ~BaseCompositeMessage, 198
  - ConstructPDV, 198
- gdcmm::network::BasePDU, 198
  - ~BasePDU, 200
  - IsLastFragment, 200
  - Print, 200
  - Read, 200
  - Size, 200
  - Write, 200
- gdcmm::network::CEchoRQ, 233
  - AffectedSOPClassUID, 234
  - ConstructPDV, 234
  - MessageID, 234
- gdcmm::network::CEchoRSP, 234
  - ConstructPDVByDataSet, 235
- gdcmm::network::CFind, 235
- gdcmm::network::CFindCancelRQ, 236
  - ConstructPDVByDataSet, 236
- gdcmm::network::CFindRQ, 237
  - ConstructPDV, 238
- gdcmm::network::CFindRSP, 238
  - ConstructPDVByDataSet, 239
- gdcmm::network::CMoveCancelRq, 239
  - ConstructPDVByDataSet, 240
- gdcmm::network::CMoveRQ, 241
  - ConstructPDV, 241
- gdcmm::network::CMoveRSP, 242
  - ConstructPDVByDataSet, 243
- gdcmm::network::CStoreRQ, 277
  - ConstructPDV, 278
- gdcmm::network::CStoreRSP, 278
  - ConstructPDV, 279
- gdcmm::network::CompositeMessageFactory, 252
  - ConstructCEchoRQ, 253
  - ConstructCFindRQ, 253
  - ConstructCMoveRQ, 253
  - ConstructCStoreRQ, 253

- ConstructCStoreRSP, [253](#)
- gdcmm::network::DIMSE, [325](#)
  - C\_CANCEL\_RQ, [326](#)
  - C\_ECHO\_RQ, [326](#)
  - C\_ECHO\_RSP, [326](#)
  - C\_FIND\_RQ, [325](#)
  - C\_FIND\_RSP, [325](#)
  - C\_GET\_RQ, [325](#)
  - C\_GET\_RSP, [325](#)
  - C\_MOVE\_RQ, [325](#)
  - C\_MOVE\_RSP, [326](#)
  - C\_STORE\_RQ, [325](#)
  - C\_STORE\_RSP, [325](#)
  - CommandTypes, [325](#)
  - N\_ACTION\_RQ, [326](#)
  - N\_ACTION\_RSP, [326](#)
  - N\_CREATE\_RQ, [326](#)
  - N\_CREATE\_RSP, [326](#)
  - N\_DELETE\_RQ, [326](#)
  - N\_DELETE\_RSP, [326](#)
  - N\_EVENT\_REPORT\_RQ, [326](#)
  - N\_EVENT\_REPORT\_RSP, [326](#)
  - N\_GET\_RQ, [326](#)
  - N\_GET\_RSP, [326](#)
  - N\_SET\_RQ, [326](#)
  - N\_SET\_RSP, [326](#)
- gdcmm::network::ImplementationClassUIDSub, [454](#)
  - ImplementationClassUIDSub, [454](#)
  - Print, [454](#)
  - Read, [454](#)
  - Size, [454](#)
  - Write, [454](#)
- gdcmm::network::ImplementationUIDSub, [454](#)
  - ImplementationUIDSub, [455](#)
  - Write, [455](#)
- gdcmm::network::ImplementationVersionNameSub, [455](#)
  - ImplementationVersionNameSub, [455](#)
  - Print, [455](#)
  - Read, [455](#)
  - Size, [455](#)
  - Write, [455](#)
- gdcmm::network::MaximumLengthSub, [503](#)
  - GetMaximumLength, [504](#)
  - MaximumLengthSub, [504](#)
  - Print, [504](#)
  - Read, [504](#)
  - SetMaximumLength, [504](#)
  - Size, [504](#)
  - Write, [504](#)
- gdcmm::network::PDUFactory, [564](#)
  - ConstructAbortPDU, [565](#)
  - ConstructPDU, [565](#)
  - ConstructReleasePDU, [565](#)
  - CreateCEchoPDU, [565](#)
  - CreateCFindPDU, [565](#)
  - CreateCMovePDU, [565](#)
  - CreateCStoreRQPDU, [565](#)
  - CreateCStoreRSPPDU, [565](#)
  - DetermineEventByPDU, [565](#)
  - GetPDVs, [565](#)
- gdcmm::network::PDataTFPDU, [556](#)
  - AddPresentationDataValue, [558](#)
  - GetNumberOfPresentationDataValues, [558](#)
  - GetPresentationDataValue, [558](#)
  - IsLastFragment, [558](#)
  - PDataTFPDU, [558](#)
  - Print, [558](#)
  - Read, [558](#)
  - ReadInto, [558](#)
  - Size, [558](#)
  - SizeType, [558](#)
  - Write, [558](#)
- gdcmm::network::PresentationContextAC, [594](#)
  - GetPresentationContextID, [594](#)
  - GetReason, [595](#)
  - GetTransferSyntax, [595](#)
  - PresentationContextAC, [594](#)
  - Print, [595](#)
  - Read, [595](#)
  - SetPresentationContextID, [595](#)
  - SetReason, [595](#)
  - SetTransferSyntax, [595](#)
  - Size, [595](#)
  - Write, [595](#)
- gdcmm::network::PresentationContextRQ, [597](#)
  - AddTransferSyntax, [598](#)
  - GetAbstractSyntax, [598](#)
  - GetNumberOfTransferSyntaxes, [598](#)
  - GetPresentationContextID, [598](#)
  - GetTransferSyntax, [598, 599](#)
  - GetTransferSyntaxes, [599](#)
  - operator==, [599](#)
  - PresentationContextRQ, [598](#)
  - Print, [599](#)
  - Read, [599](#)
  - SetAbstractSyntax, [599](#)
  - SetPresentationContextID, [599](#)
  - Size, [599](#)
  - SizeType, [598](#)
  - Write, [599](#)
- gdcmm::network::PresentationDataValue, [599](#)
  - ConcatenatePDVBlobs, [600](#)
  - ConcatenatePDVBlobsAsExplicit, [600](#)
  - GetBlob, [600](#)
  - GetIsCommand, [600](#)
  - GetIsLastFragment, [600](#)
  - GetMessageHeader, [600](#)
  - GetPresentationContextID, [600](#)

- PresentationDataValue, 600
- Print, 600
- Read, 600
- ReadInto, 600
- SetBlob, 600
- SetCommand, 600
- SetDataSet, 600
- SetLastFragment, 601
- SetMessageHeader, 601
- SetPresentationContextID, 601
- Size, 601
- Write, 601
- gdcmm::network::RoleSelectionSub, 640
  - Print, 641
  - Read, 641
  - RoleSelectionSub, 641
  - SetTuple, 641
  - Size, 641
  - Write, 641
- gdcmm::network::SOPClassExtendedNegociationSub, 688
  - Print, 689
  - Read, 689
  - SOPClassExtendedNegociationSub, 689
  - SetTuple, 689
  - Size, 689
  - Write, 689
- gdcmm::network::ServiceClassApplicationInformation, 673
  - Print, 673
  - Read, 673
  - ServiceClassApplicationInformation, 673
  - SetTuple, 673
  - Size, 674
  - Write, 674
- gdcmm::network::TableRow, 740
  - ~TableRow, 741
  - TableRow, 741
  - transitions, 741
- gdcmm::network::TransferSyntaxSub, 761
  - GetName, 762
  - operator==, 762
  - Print, 762
  - Read, 762
  - SetName, 762
  - SetNameFromUID, 762
  - Size, 762
  - TransferSyntaxSub, 762
  - Write, 762
- gdcmm::network::Transition, 762
  - ~Transition, 763
  - mAction, 763
  - mEnd, 763
  - MakeNew, 763
  - Transition, 763
- gdcmm::network::ULAction, 786
  - ~ULAction, 788
  - PerformAction, 788
  - ULAction, 788
- gdcmm::network::ULActionAA1, 789
  - PerformAction, 789
- gdcmm::network::ULActionAA2, 790
  - PerformAction, 790
- gdcmm::network::ULActionAA3, 791
  - PerformAction, 791
- gdcmm::network::ULActionAA4, 792
  - PerformAction, 792
- gdcmm::network::ULActionAA5, 793
  - PerformAction, 793
- gdcmm::network::ULActionAA6, 794
  - PerformAction, 794
- gdcmm::network::ULActionAA7, 795
  - PerformAction, 795
- gdcmm::network::ULActionAA8, 796
  - PerformAction, 796
- gdcmm::network::ULActionAE1, 797
  - PerformAction, 797
- gdcmm::network::ULActionAE2, 798
  - PerformAction, 798
- gdcmm::network::ULActionAE3, 799
  - PerformAction, 799
- gdcmm::network::ULActionAE4, 800
  - PerformAction, 800
- gdcmm::network::ULActionAE5, 801
  - PerformAction, 801
- gdcmm::network::ULActionAE6, 802
  - PerformAction, 802
- gdcmm::network::ULActionAE7, 803
  - PerformAction, 803
- gdcmm::network::ULActionAE8, 804
  - PerformAction, 804
- gdcmm::network::ULActionAR1, 805
  - PerformAction, 805
- gdcmm::network::ULActionAR10, 806
  - PerformAction, 806
- gdcmm::network::ULActionAR2, 807
  - PerformAction, 807
- gdcmm::network::ULActionAR3, 808
  - PerformAction, 808
- gdcmm::network::ULActionAR4, 809
  - PerformAction, 809
- gdcmm::network::ULActionAR5, 810
  - PerformAction, 810
- gdcmm::network::ULActionAR6, 811
  - PerformAction, 811
- gdcmm::network::ULActionAR7, 812
  - PerformAction, 812
- gdcmm::network::ULActionAR8, 813
  - PerformAction, 813
- gdcmm::network::ULActionAR9, 814
  - PerformAction, 814



- PerformAction, [814](#)
- gdcmm::network::ULActionDT1, [815](#)
  - PerformAction, [815](#)
- gdcmm::network::ULActionDT2, [816](#)
  - PerformAction, [816](#)
- gdcmm::network::ULBasicCallback, [817](#)
  - ~ULBasicCallback, [818](#)
  - GetDataSets, [818](#)
  - GetResponses, [818](#)
  - HandleDataSet, [818](#)
  - HandleResponse, [818](#)
  - ULBasicCallback, [818](#)
- gdcmm::network::ULConnection, [818](#)
  - ~ULConnection, [819](#)
  - AddAcceptedPresentationContext, [819](#)
  - FindContext, [819](#)
  - GetAcceptedPresentationContexts, [820](#)
  - GetConnectionInfo, [820](#)
  - GetMaxPDUSize, [820](#)
  - GetPresentationContextACByID, [820](#)
  - GetPresentationContextIDFromPresentationContext, [820](#)
  - GetPresentationContextRQByID, [820](#)
  - GetPresentationContexts, [820](#)
  - GetProtocol, [820](#)
  - GetState, [820](#)
  - GetTimer, [820](#)
  - InitializeConnection, [820](#)
  - InitializeIncomingConnection, [820](#)
  - SetMaxPDUSize, [820](#)
  - SetPresentationContexts, [820](#)
  - SetState, [820](#)
  - StopProtocol, [820](#)
  - ULActionAE6, [820](#)
  - ULConnection, [819](#)
  - ULConnectionManager, [821](#)
- gdcmm::network::ULConnectionCallback, [821](#)
  - ~ULConnectionCallback, [822](#)
  - DataSetHandled, [822](#)
  - DataSetHandles, [822](#)
  - HandleDataSet, [822](#)
  - HandleResponse, [822](#)
  - mImplicit, [822](#)
  - ResetHandledDataSet, [822](#)
  - SetImplicitFlag, [822](#)
  - ULConnectionCallback, [822](#)
- gdcmm::network::ULConnectionInfo, [822](#)
  - GetCalledAETitle, [823](#)
  - GetCalledComputerName, [823](#)
  - GetCalledIPAddress, [823](#)
  - GetCalledIPPort, [823](#)
  - GetCallingAETitle, [823](#)
  - GetMaxPDULength, [823](#)
  - Initialize, [823](#)
  - SetMaxPDULength, [823](#)
  - ULConnectionInfo, [823](#)
- gdcmm::network::ULConnectionManager, [824](#)
  - ~ULConnectionManager, [825](#)
  - BreakConnection, [825](#)
  - BreakConnectionNow, [825](#)
  - EstablishConnection, [825](#)
  - EstablishConnectionMove, [825](#)
  - SendEcho, [826](#)
  - SendFind, [826](#)
  - SendMove, [826](#)
  - SendStore, [826](#)
  - ULConnectionManager, [825](#)
- gdcmm::network::ULEvent, [826](#)
  - ~ULEvent, [827](#)
  - GetEvent, [827](#)
  - GetPDUs, [827](#)
  - SetEvent, [827](#)
  - SetPDU, [827](#)
  - ULEvent, [827](#)
- gdcmm::network::ULTransitionTable, [827](#)
  - HandleEvent, [828](#)
  - PrintTable, [828](#)
  - ULTransitionTable, [828](#)
- gdcmm::network::ULWritingCallback, [828](#)
  - ~ULWritingCallback, [829](#)
  - HandleDataSet, [829](#)
  - HandleResponse, [829](#)
  - SetDirectory, [829](#)
  - ULWritingCallback, [829](#)
- gdcmm::network::UserInformation, [837](#)
  - ~UserInformation, [837](#)
  - AddRoleSelectionSub, [837](#)
  - AddSOPClassExtendedNegotiationSub, [837](#)
  - GetMaximumLengthSub, [837](#)
  - operator=, [838](#)
  - Print, [838](#)
  - Read, [838](#)
  - Size, [838](#)
  - UserInformation, [837](#)
  - Write, [838](#)
- gdcmm::static\_assert\_test< x >, [699](#)
- gdcmm::terminal, [136](#)
  - Attribute, [137](#)
  - black, [137](#)
  - blink, [137](#)
  - blue, [138](#)
  - bright, [137](#)
  - CONSOLE, [138](#)
  - Color, [137](#)
  - cyan, [138](#)
  - dim, [137](#)
  - green, [137](#)
  - hidden, [137](#)



- magenta, [138](#)
- Mode, [138](#)
- red, [137](#)
- reset, [137](#)
- reverse, [137](#)
- setAttribute, [138](#)
- setbgcolor, [138](#)
- setfgcolor, [138](#)
- setmode, [138](#)
- underline, [137](#)
- VT100, [138](#)
- white, [138](#)
- yellow, [137](#)
- gdcmAAAbortPDU.h, [929](#)
- gdcmAAAssociateACPDU.h, [930](#)
- gdcmAAAssociateRJPDU.h, [930](#)
- gdcmAAAssociateRQPDU.h, [931](#)
- gdcmARTIMTimer.h, [938](#)
- gdcmAReleaseRPPDU.h, [936](#)
- gdcmAReleaseRQPDU.h, [937](#)
- gdcmASN1.h, [939](#)
- gdcmAbstractSyntax.h, [932](#)
- gdcmAnonymizeEvent.h, [933](#)
- gdcmAnonymizer.h, [934](#)
- gdcmApplicationContext.h, [935](#)
- gdcmApplicationEntity.h, [936](#)
- gdcmAssertAlwaysMacro
  - gdcmTrace.h, [1155](#)
- gdcmAssertMacro
  - gdcmTrace.h, [1155](#)
- gdcmAsynchronousOperationsWindowSub.h, [939](#)
- gdcmAttribute.h, [940](#)
- gdcmAudioCodec.h, [942](#)
- gdcmBase64.h, [942](#)
- gdcmBaseCompositeMessage.h, [943](#)
- gdcmBasePDU.h, [944](#)
- gdcmBaseRootQuery.h, [945](#)
- gdcmBasicOffsetTable.h, [946](#)
- gdcmBitmap.h, [948](#)
- gdcmBitmapToBitmapFilter.h, [949](#)
- gdcmBoxRegion.h, [949](#)
- gdcmByteBuffer.h, [950](#)
- gdcmByteSwap.h, [951](#)
- gdcmByteSwapFilter.h, [952](#)
- gdcmByteValue.h, [953](#)
- gdcmCAPICryptoFactory.h, [954](#)
- gdcmCAPICryptographicMessageSyntax.h, [954](#)
- gdcmCEchoMessages.h, [955](#)
- gdcmCFindMessages.h, [956](#)
- gdcmCMoveMessages.h, [957](#)
- gdcmCP246ExplicitDataElement.h, [965](#)
- gdcmCSAElement.h, [967](#)
- gdcmCSAHeader.h, [968](#)
- gdcmCSAHeaderDictEntry.h, [969](#)
- gdcmCSAHeaderDictEntry.h, [971](#)
- gdcmCStoreMessages.h, [972](#)
- gdcmCodeString.h, [959](#)
- gdcmCodec.h, [957](#)
- gdcmCoder.h, [958](#)
- gdcmCommand.h, [960](#)
- gdcmCommandDataSet.h, [962](#)
- gdcmCompositeMessageFactory.h, [962](#)
- gdcmCompositeNetworkFunctions.h, [963](#)
- gdcmConstCharWrapper.h, [964](#)
- gdcmCryptoFactory.h, [965](#)
- gdcmCryptographicMessageSyntax.h, [966](#)
- gdcmCurve.h, [973](#)
- gdcmDICOMDIR.h, [983](#)
- gdcmDICOMDIRGenerator.h, [984](#)
- gdcmDIMSE.h, [991](#)
- gdcmDataElement.h, [974](#)
- gdcmDataEvent.h, [975](#)
- gdcmDataSet.h, [976](#)
- gdcmDataSetEvent.h, [977](#)
- gdcmDataSetHelper.h, [978](#)
- gdcmDebugMacro
  - gdcmTrace.h, [1155](#)
- gdcmDecoder.h, [979](#)
- gdcmDefinedTerms.h, [980](#)
- gdcmDeflateStream.h, [981](#)
- gdcmDefs.h, [981](#)
- gdcmDeltaEncodingCodec.h, [983](#)
- gdcmDict.h, [985](#)
- gdcmDictConverter.h, [987](#)
- gdcmDictEntry.h, [987](#)
- gdcmDictPrinter.h, [989](#)
- gdcmDicts.h, [989](#)
- gdcmDirectionCosines.h, [991](#)
- gdcmDirectory.h, [992](#)
- gdcmDirectoryHelper.h, [993](#)
- gdcmDummyValueGenerator.h, [994](#)
- gdcmDumper.h, [994](#)
- gdcmElement.h, [995](#)
- gdcmEncapsulatedDocument.h, [997](#)
- gdcmEnumeratedValues.h, [997](#)
- gdcmErrorMacro
  - gdcmTrace.h, [1156](#)
- gdcmEvent.h, [998](#)
- gdcmEventMacro, [999](#)
- gdcmEventMacro
  - gdcmEvent.h, [999](#)
- gdcmException.h, [1000](#)
- gdcmExplicitDataElement.h, [1000](#)
- gdcmExplicitImplicitDataElement.h, [1001](#)
- gdcmFiducials.h, [1002](#)
- gdcmFile.h, [1003](#)
- gdcmFileAnonymizer.h, [1004](#)
- gdcmFileChangeTransferSyntax.h, [1004](#)

gdcmlFileDerivation.h, 1005  
 gdcmlFileExplicitFilter.h, 1006  
 gdcmlFileMetaInformation.h, 1007  
 gdcmlFileNameEvent.h, 1009  
 gdcmlFileSet.h, 1010  
 gdcmlFileStreamer.h, 1012  
 gdcmlFilename.h, 1008  
 gdcmlFilenameGenerator.h, 1010  
 gdcmlFindPatientRootQuery.h, 1012  
 gdcmlFindStudyRootQuery.h, 1013  
 gdcmlFragment.h, 1014  
 gdcmlGlobal.h, 1016  
 gdcmlGroupDict.h, 1017  
 gdcmlIOD.h, 1035  
 gdcmlIODEntry.h, 1036  
 gdcmlIODs.h, 1038  
 gdcmlIPPSorter.h, 1039  
 gdcmlIconImage.h, 1017  
 gdcmlIconImageFilter.h, 1018  
 gdcmlIconImageGenerator.h, 1019  
 gdcmlImage.h, 1020  
 gdcmlImageApplyLookupTable.h, 1021  
 gdcmlImageChangePhotometricInterpretation.h, 1021  
 gdcmlImageChangePlanarConfiguration.h, 1022  
 gdcmlImageChangeTransferSyntax.h, 1023  
 gdcmlImageCodec.h, 1024  
 gdcmlImageConverter.h, 1025  
 gdcmlImageFragmentSplitter.h, 1026  
 gdcmlImageHelper.h, 1027  
 gdcmlImageReader.h, 1028  
 gdcmlImageRegionReader.h, 1028  
 gdcmlImageToImageFilter.h, 1029  
 gdcmlImageWriter.h, 1030  
 gdcmlImplementationClassUIDSub.h, 1031  
 gdcmlImplementationUIDSub.h, 1032  
 gdcmlImplementationVersionNameSub.h, 1033  
 gdcmlImplicitDataElement.h, 1034  
 gdcmlItem.h, 1040  
 gdcmlJPEG12Codec.h, 1042  
 gdcmlJPEG16Codec.h, 1042  
 gdcmlJPEG2000Codec.h, 1043  
 gdcmlJPEG8Codec.h, 1044  
 gdcmlJPEGCodec.h, 1045  
 gdcmlJPEGLSCodec.h, 1046  
 gdcmlJSON.h, 1047  
 gdcmlKAKADUCodec.h, 1048  
 gdcmlLO.h, 1049  
 gdcmlLegacyMacro.h, 1048  
     GDCM\_LEGACY, 1049  
     GDCM\_LEGACY\_BODY, 1049  
     GDCM\_LEGACY\_REPLACED\_BODY, 1049  
 gdcmlLookupTable.h, 1050  
 gdcmlMD5.h, 1057  
 gdcmlMacro.h, 1051  
 gdcmlMacroEntry.h, 1053  
     GDCMMACROENTRY\_H, 1054  
 gdcmlMacros.h, 1054  
 gdcmlMaximumLengthSub.h, 1056  
 gdcmlMediaStorage.h, 1058  
 gdcmlMeshPrimitive.h, 1059  
 gdcmlModule.h, 1060  
 gdcmlModuleEntry.h, 1062  
 gdcmlModules.h, 1063  
 gdcmlMovePatientRootQuery.h, 1065  
 gdcmlMoveStudyRootQuery.h, 1066  
 gdcmlNestedModuleEntries.h, 1066  
 gdcmlNetworkEvents.h, 1068  
 gdcmlNetworkStateID.h, 1069  
 gdcmlObject.h, 1070  
 gdcmlOpenSSLCryptoFactory.h, 1070  
 gdcmlOpenSSLCryptographicMessageSyntax.h, 1071  
 gdcmlOpenSSL7CryptoFactory.h, 1072  
 gdcmlOpenSSL7CryptographicMessageSyntax.h, 1073  
 gdcmlOrientation.h, 1075  
 gdcmlOverlay.h, 1075  
 gdcmlPDBelement.h, 1080  
 gdcmlPDBHeader.h, 1082  
 gdcmlPDFCodec.h, 1082  
 gdcmlPDUFactory.h, 1083  
 gdcmlPDataTFPDU.h, 1079  
 gdcmlPGXCodec.h, 1084  
 gdcmlPNMCodec.h, 1091  
 gdcmlPVRGCodec.h, 1100  
 gdcmlParseException.h, 1076  
 gdcmlParser.h, 1078  
 gdcmlPatient.h, 1078  
 gdcmlPersonName.h, 1084  
 gdcmlPhotometricInterpretation.h, 1085  
 gdcmlPixelFormat.h, 1086  
 gdcmlPixmap.h, 1087  
 gdcmlPixmapReader.h, 1088  
 gdcmlPixmapToPixmapFilter.h, 1089  
 gdcmlPixmapWriter.h, 1090  
 gdcmlPreamble.h, 1092  
 gdcmlPresentationContext.h, 1093  
 gdcmlPresentationContextAC.h, 1094  
 gdcmlPresentationContextGenerator.h, 1095  
 gdcmlPresentationContextRQ.h, 1096  
 gdcmlPresentationDataValue.h, 1097  
 gdcmlPrinter.h, 1097  
 gdcmlPrivateTag.h, 1098  
 gdcmlProgressEvent.h, 1100  
 gdcmlPythonFilter.h, 1101  
 gdcmlQueryBase.h, 1102  
 gdcmlQueryFactory.h, 1103  
 gdcmlQueryImage.h, 1104  
 gdcmlQueryPatient.h, 1105  
 gdcmlQuerySeries.h, 1106

gdcmQueryStudy.h, 1107  
gdcmRAWCodec.h, 1107  
gdcmRLECodec.h, 1111  
gdcmReader.h, 1108  
gdcmRegion.h, 1110  
gdcmRescaler.h, 1111  
gdcmRoleSelectionSub.h, 1112  
gdcmSHA1.h, 1125  
gdcmSOPClassExtendedNegociationSub.h, 1128  
gdcmSOPClassUIDToIOD.h, 1129  
gdcmScanner.h, 1113  
gdcmSegment.h, 1114  
gdcmSegmentHelper.h, 1116  
gdcmSegmentReader.h, 1117  
gdcmSegmentWriter.h, 1118  
gdcmSegmentedPaletteColorLookupTable.h, 1115  
gdcmSequenceOfFragments.h, 1119  
gdcmSequenceOfItems.h, 1120  
gdcmSerieHelper.h, 1121  
gdcmSeries.h, 1122  
gdcmServiceClassApplicationInformation.h, 1124  
gdcmServiceClassUser.h, 1125  
gdcmSimpleSubjectWatcher.h, 1126  
gdcmSmartPointer.h, 1127  
gdcmSorter.h, 1130  
gdcmSpacing.h, 1132  
gdcmSpectroscopy.h, 1132  
gdcmSplitMosaicFilter.h, 1133  
gdcmStaticAssert.h, 1134  
    GDCM\_DO\_JOIN, 1134  
    GDCM\_DO\_JOIN2, 1134  
    GDCM\_JOIN, 1134  
    GDCM\_STATIC\_ASSERT, 1134  
gdcmStreamImageReader.h, 1135  
gdcmStreamImageWriter.h, 1135  
gdcmString.h, 1136  
gdcmStringFilter.h, 1137  
gdcmStudy.h, 1138  
gdcmSubject.h, 1139  
gdcmSurface.h, 1139  
gdcmSurfaceHelper.h, 1141  
gdcmSurfaceReader.h, 1141  
gdcmSurfaceWriter.h, 1142  
gdcmSwapCode.h, 1143  
gdcmSwapper.h, 1144  
gdcmSystem.h, 1145  
gdcmTable.h, 1146  
gdcmTableEntry.h, 1146  
gdcmTableReader.h, 1147  
gdcmTag.h, 1149  
gdcmTagPath.h, 1149  
gdcmTagToVR.h, 1150  
gdcmTerminal.h, 1151  
gdcmTestDriver.h, 1153  
gdcmTesting.h, 1153  
gdcmTrace.h, 1154  
    GDCM\_FUNCTION, 1155  
    gdcmAssertAlwaysMacro, 1155  
    gdcmAssertMacro, 1155  
    gdcmDebugMacro, 1155  
    gdcmErrorMacro, 1156  
    gdcmWarningMacro, 1156  
gdcmTransferSyntax.h, 1157  
gdcmTransferSyntaxSub.h, 1158  
gdcmType.h, 1159  
gdcmTypes.h, 1160  
gdcmUIDGenerator.h, 1161  
gdcmUIDs.h, 1161  
gdcmULAction.h, 1162  
gdcmULActionAA.h, 1163  
gdcmULActionAE.h, 1164  
gdcmULActionAR.h, 1165  
gdcmULActionDT.h, 1166  
gdcmULBasicCallback.h, 1166  
gdcmULConnection.h, 1167  
gdcmULConnectionCallback.h, 1168  
gdcmULConnectionInfo.h, 1169  
gdcmULConnectionManager.h, 1171  
gdcmULEvent.h, 1171  
gdcmULTransitionTable.h, 1172  
gdcmULWritingCallback.h, 1174  
gdcmUNExplicitDataElement.h, 1174  
gdcmUNExplicitImplicitDataElement.h, 1175  
gdcmUUIIDGenerator.h, 1179  
gdcmUnpacker12Bits.h, 1175  
gdcmUsage.h, 1176  
gdcmUserInformation.h, 1178  
gdcmVL.h, 1183  
gdcmVM.h, 1184  
    TYPETOLENGTH, 1185  
gdcmVR.h, 1185  
    TYPETOENCODING, 1187  
    VRTypeTemplateCase, 1187  
gdcmVR16ExplicitDataElement.h, 1187  
gdcmValidate.h, 1179  
gdcmValue.h, 1180  
gdcmValueIO.h, 1181  
gdcmVersion.h, 1182  
gdcmWarningMacro  
    gdcmTrace.h, 1156  
gdcmWaveform.h, 1188  
gdcmWin32.h, 1188  
    GDCM\_EXPORT, 1188  
gdcmWriter.h, 1189  
gdcmXMLDictReader.h, 1190  
gdcmXMLPrinter.h, 1190  
gdcmXMLPrivateDictReader.h, 1191  
gdcmanon.dox, 933

- gdcmconv.dox, [964](#)
- gdcmdiff.dox, [990](#)
- gdcmdump.dox, [994](#)
- gdcmgendir.dox, [1016](#)
- gdcmimg.dox, [1031](#)
- gdcminfo.dox, [1034](#)
- gdcmpap3.dox, [1076](#)
- gdcmpdf.dox, [1082](#)
- gdcmraw.dox, [1107](#)
- gdcmscanner.dox, [1113](#)
- gdcm SCU.dox, [1114](#)
- gdcm tar.dox, [1151](#)
- gdcmviewer.dox, [1183](#)
- gdcmxml.dox, [1190](#)
- GeneralECGWaveformStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [775](#)
- GeneralElectricMagneticResonanceImageStorage
  - gdcm::MediaStorage, [509](#)
- GeneralPurposePerformedProcedureStepSOPClass
  - gdcm::UIDs, [777](#)
- GeneralPurposeScheduledProcedureStepSOPClass
  - gdcm::UIDs, [777](#)
- GeneralPurposeWorklistInformationModelFIND
  - gdcm::UIDs, [777](#)
- GeneralPurposeWorklistManagementMetaSOPClass
  - gdcm::UIDs, [777](#)
- GeneralRelevantPatientInformationQuery
  - gdcm::UIDs, [777](#)
- Generate
  - gdcm::DICOMDIRGenerator, [312](#)
  - gdcm::DummyValueGenerator, [332](#)
  - gdcm::FilenameGenerator, [390](#)
  - gdcm::IconImageGenerator, [411](#)
  - gdcm::UIDGenerator, [767](#)
  - gdcm::UUIDGenerator, [838](#)
- GenerateFromFilenames
  - gdcm::PresentationContextGenerator, [596](#)
- GenerateFromUID
  - gdcm::PresentationContextGenerator, [597](#)
- GenerateUUID
  - gdcm::UIDGenerator, [767](#)
- Get
  - gdcm::ByteBuffer, [222](#)
- GetAETitle
  - gdcm::ServiceClassUser, [676](#)
- GetALGOType
  - gdcm::Segment, [650](#)
- GetALGOTypeString
  - gdcm::Segment, [650](#)
- GetAbbreviation
  - gdcm::GroupDict, [407](#)
- GetAbstractSyntax
  - gdcm::PresentationContext, [593](#)
  - gdcm::network::PresentationContextRQ, [598](#)
- GetAbstractSyntaxUID
  - gdcm::BaseRootQuery, [203](#)
  - gdcm::FindPatientRootQuery, [398](#)
  - gdcm::FindStudyRootQuery, [401](#)
  - gdcm::MovePatientRootQuery, [528](#)
  - gdcm::MoveStudyRootQuery, [530](#)
- GetAcceptedPresentationContexts
  - gdcm::network::ULConnection, [820](#)
- GetAlgorithmFamily
  - gdcm::Surface, [719](#)
- GetAlgorithmName
  - gdcm::Surface, [719](#)
- GetAlgorithmVersion
  - gdcm::Surface, [719](#)
- GetAllFilenamesFromTagToValue
  - gdcm::Scanner, [646](#)
- GetAllRequiredTags
  - gdcm::QueryBase, [614](#)
- GetAllTags
  - gdcm::QueryBase, [614](#)
- GetAnatomicRegion
  - gdcm::Segment, [650](#)
- GetAsDataElement
  - gdcm::Attribute, [171](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [184](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [340](#)
  - gdcm::PrivateTag, [607](#)
  - gdcm::network::AbstractSyntax, [151](#)
- GetAsPoints
  - gdcm::Curve, [281](#)
- GetAsString
  - gdcm::CodeString, [247](#)
- GetAxisOfRotation
  - gdcm::Surface, [719](#)
- GetBasicApplicationLevelConfidentialityProfileAttributes
  - gdcm::Anonymizer, [157](#)
- GetBitPosition
  - gdcm::Overlay, [549](#)
- GetBitSample
  - gdcm::LookupTable, [499](#)
- GetBitsAllocated
  - gdcm::Overlay, [549](#)
  - gdcm::PixelFormat, [574](#)
- GetBitsStored
  - gdcm::PixelFormat, [574](#)
- GetBlob
  - gdcm::network::PresentationDataValue, [600](#)
- GetBuffer
  - gdcm::Bitmap, [212](#)

- gdcm::ByteValue, [227](#)
- gdcm::Parser, [555](#)
- gdcm::SequenceOfFragments, [662](#)
- GetBuffer2
  - gdcm::Bitmap, [212](#)
- GetBufferAsRGBA
  - gdcm::LookupTable, [499](#)
- GetBufferLength
  - gdcm::Bitmap, [212](#)
  - gdcm::JPEGLSCodec, [489](#)
  - gdcm::PNMCodec, [590](#)
  - gdcm::RLECodec, [639](#)
- GetBuildVersion
  - gdcm::Version, [843](#)
- GetByteValue
  - gdcm::CSAElement, [266](#)
  - gdcm::DataElement, [286](#)
- GetCSADataInfo
  - gdcm::CSAHeader, [271](#)
- GetCSAEEnd
  - gdcm::CSAHeader, [271](#)
- GetCSAElementByName
  - gdcm::CSAHeader, [271](#)
- GetCSAHeaderDict
  - gdcm::Dicts, [324](#)
- GetCSAHeaderDictEntry
  - gdcm::CSAHeaderDict, [274](#)
- GetCSAImageHeaderInfoTag
  - gdcm::CSAHeader, [271](#)
- GetCSASeriesHeaderInfoTag
  - gdcm::CSAHeader, [271](#)
- GetCTImageSeriesUIDs
  - gdcm::DirectoryHelper, [331](#)
- GetCWD
  - gdcm::System, [735](#)
- GetCalledAETitle
  - gdcm::ServiceClassUser, [676](#)
  - gdcm::network::AAssociateRQPDU, [147](#)
  - gdcm::network::ULConnectionInfo, [823](#)
- GetCalledComputerName
  - gdcm::network::ULConnectionInfo, [823](#)
- GetCalledIPAddress
  - gdcm::network::ULConnectionInfo, [823](#)
- GetCalledIPPort
  - gdcm::network::ULConnectionInfo, [823](#)
- GetCallingAETitle
  - gdcm::network::AAssociateRQPDU, [148](#)
  - gdcm::network::ULConnectionInfo, [823](#)
- GetCenterOfRotation
  - gdcm::Surface, [719](#)
- GetCharacterFromCurrentLocale
  - gdcm::QueryFactory, [616](#)
- GetCipherType
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [540](#)
  - gdcm::OpenSSLP7CryptographicMessageSyntax, [543](#)
- GetCodec
  - gdcm::FileChangeTransferSyntax, [374](#)
- GetColorLevel
  - vtkImageColorViewer, [895](#)
- GetColorWindow
  - vtkImageColorViewer, [895](#)
- GetColumns
  - gdcm::Bitmap, [213](#)
  - gdcm::Overlay, [549](#)
- GetCommand
  - gdcm::Subject, [715](#)
- GetConnectionInfo
  - gdcm::network::ULConnection, [820](#)
- GetConstructorString
  - gdcm::Dicts, [323](#)
- GetContourReferencedFrameOfReferenceClassUID
  - vtkRTStructSetProperties, [914](#)
- GetContourReferencedFrameOfReferenceInstanceUID
  - vtkRTStructSetProperties, [914](#)
- GetCryptographicMessageSyntax
  - gdcm::Anonymizer, [157](#)
- GetCurrentByteIndex
  - gdcm::Parser, [555](#)
- GetCurrentDateTime
  - gdcm::System, [734](#)
- GetCurrentModuleFileName
  - gdcm::System, [734](#)
- GetCurrentProcessFileName
  - gdcm::System, [735](#)
- GetCurrentResourcesDirectory
  - gdcm::System, [735](#)
- GetCurve
  - gdcm::Pixmap, [578](#), [579](#)
- GetCurveDataDescriptor
  - gdcm::Curve, [281](#)
- GetDEEnd
  - gdcm::DataSet, [299](#)
- GetDES
  - gdcm::DataSet, [299](#)
- GetData
  - gdcm::DataEvent, [294](#)
- GetDataElement
  - gdcm::Bitmap, [213](#)
  - gdcm::DataSet, [298](#), [299](#)
  - gdcm::Item, [469](#)
- GetDataExtraRoot
  - gdcm::Testing, [751](#)
- GetDataLength
  - gdcm::DataEvent, [294](#)
- GetDataRoot

- gdcmm::Testing, 751
- GetDataSet
  - gdcmm::CSAHeader, 271
  - gdcmm::DataSetEvent, 303
  - gdcmm::File, 367, 368
- GetDataSetTransferSyntax
  - gdcmm::FileMetaInformation, 381
- GetDataSets
  - gdcmm::network::ULBasicCallback, 818
- GetDataValueRepresentation
  - gdcmm::Curve, 282
- GetDebugFlag
  - gdcmm::Trace, 755
- GetDebugStream
  - gdcmm::Trace, 755
- GetDecodeLength
  - gdcmm::Base64, 196
- GetDefaultTransferSyntax
  - gdcmm::PresentationContextGenerator, 597
- GetDefs
  - gdcmm::Global, 405
  - gdcmm::TableReader, 739
- GetDescription
  - gdcmm::CSAHeaderDictEntry, 275
  - gdcmm::Exception, 360
  - gdcmm::ModuleEntry, 524
  - gdcmm::Overlay, 549
- GetDescriptiveName
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
  - vtkGDCMImageWriter, 873
- GetDict
  - gdcmm::XMLDictReader, 923
- GetDictEntry
  - gdcmm::Dict, 315
  - gdcmm::Dicts, 324
  - gdcmm::PrivateDict, 605
- GetDictEntryByKeyword
  - gdcmm::Dict, 315
- GetDictEntryByName
  - gdcmm::Dict, 315
- GetDictName
  - gdcmm::DictConverter, 317
- GetDictVM
  - gdcmm::Attribute, 171
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 178
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, 184
- GetDictVR
  - gdcmm::Attribute, 171
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 178
- gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, 185
- GetDicts
  - gdcmm::Global, 405
- GetDimension
  - gdcmm::Bitmap, 213
- GetDimensions
  - gdcmm::Bitmap, 213
  - gdcmm::Curve, 282
  - gdcmm::ImageCodec, 432
- GetDimensionsValue
  - gdcmm::ImageHelper, 441
- GetDimensionsValueForResolution
  - gdcmm::StreamImageReader, 701
- GetDirectionCosines
  - gdcmm::Image, 415
- GetDirectionCosinesFromDataSet
  - gdcmm::ImageHelper, 441
- GetDirectionCosinesTolerance
  - gdcmm::IPPSorter, 465
- GetDirectionCosinesValue
  - gdcmm::ImageHelper, 441
- GetDirectories
  - gdcmm::Directory, 329
- GetElapsedTime
  - gdcmm::network::ARTIMTimer, 167
- GetElement
  - gdcmm::Tag, 744
- GetElementTag
  - gdcmm::Tag, 744
- GetEncodeLength
  - gdcmm::Base64, 196
- GetErrorCode
  - gdcmm::Parser, 555
- GetErrorFlag
  - gdcmm::Trace, 755
- GetErrorStream
  - gdcmm::Trace, 755
- GetErrorString
  - gdcmm::Parser, 555
- GetEvent
  - gdcmm::network::ULEvent, 827
- GetEventName
  - gdcmm::AnonymizeEvent, 153
  - gdcmm::DataEvent, 294
  - gdcmm::DataSetEvent, 303
  - gdcmm::Event, 358
  - gdcmm::FileNameEvent, 388
  - gdcmm::ProgressEvent, 609
- GetExtension
  - gdcmm::Filename, 385
- GetFactoryInstance
  - gdcmm::CryptoFactory, 261
- GetFile



- gdcm::Anonymizer, [157](#)
- gdcm::DICOMDIRGenerator, [312](#)
- gdcm::FileDerivation, [376](#)
- gdcm::FileExplicitFilter, [378](#)
- gdcm::IconImageFilter, [409](#)
- gdcm::PythonFilter, [613](#)
- gdcm::Reader, [630](#)
- gdcm::SplitMosaicFilter, [697](#)
- gdcm::StreamImageReader, [701](#)
- gdcm::StringFilter, [712](#)
- gdcm::Writer, [920](#)
- vtkGDCMMedicalImageProperties, [877](#)
- GetFileExtensions
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [868](#)
  - vtkGDCMImageWriter, [873](#)
- GetFileMetaInformationVersion
  - gdcm::FileMetaInformation, [382](#)
- GetFileName
  - gdcm::FileNameEvent, [388](#)
  - gdcm::Filename, [385](#)
  - gdcm::Testing, [751](#)
  - vtkGDCMImageWriter, [873](#)
  - vtkGDCMThreadedImageReader2, [890](#)
- GetFileNames
  - gdcm::Testing, [751](#)
- GetFilename
  - gdcm::FilenameGenerator, [390](#)
  - gdcm::TableReader, [739](#)
- GetFilenameFromTagToValue
  - gdcm::Scanner, [646](#)
- GetFilenames
  - gdcm::Directory, [330](#)
  - gdcm::FilenameGenerator, [390](#)
  - gdcm::Scanner, [646](#)
  - gdcm::Sorter, [693](#)
- GetFilenamesFromSeriesUIDs
  - gdcm::DirectoryHelper, [331](#)
- GetFiles
  - gdcm::FileSet, [392](#)
- GetFiniteVolume
  - gdcm::Surface, [720](#)
- GetFirstSingleSerieUIDFileSet
  - gdcm::SerieHelper, [672](#)
- GetForcePixelSpacing
  - gdcm::ImageHelper, [441](#)
- GetForceRescaleInterceptSlope
  - gdcm::ImageHelper, [441](#)
- GetFormat
  - gdcm::CSAHeader, [272](#)
- GetFragBuffer
  - gdcm::SequenceOfFragments, [662](#)
- GetFragment
  - gdcm::SequenceOfFragments, [662](#)
- GetFragmentSizeMax
  - gdcm::ImageFragmentSplitter, [439](#)
- GetFrameOfReference
  - gdcm::DirectoryHelper, [331](#)
- GetFullLength
  - gdcm::FileMetaInformation, [382](#)
- GetGDCMDataRoot
  - vtkGDCMTesting, [885](#)
- GetGDCMImplementationClassUID
  - gdcm::FileMetaInformation, [382](#)
- GetGDCMImplementationVersionName
  - gdcm::FileMetaInformation, [382](#)
- GetGDCMSourceApplicationEntityTitle
  - gdcm::FileMetaInformation, [382](#)
- GetGDCMUID
  - gdcm::UIDGenerator, [767](#)
- GetGroup
  - gdcm::Curve, [282](#)
  - gdcm::Overlay, [550](#)
  - gdcm::Tag, [744](#)
- GetHasExpired
  - gdcm::network::ARTIMTimer, [167](#)
- GetHeader
  - gdcm::File, [368](#)
- GetHeaderInfo
  - gdcm::ImageCodec, [433](#)
  - gdcm::JPEG12Codec, [473](#)
  - gdcm::JPEG16Codec, [475](#)
  - gdcm::JPEG2000Codec, [478](#)
  - gdcm::JPEG8Codec, [481](#)
  - gdcm::JPEGCodec, [485](#)
  - gdcm::JPEGLSCodec, [489](#)
  - gdcm::PGXCodec, [569](#)
  - gdcm::PNMCodec, [590](#)
  - gdcm::RAWCodec, [627](#)
  - gdcm::RLECodec, [639](#)
- GetHierarchicalSearchTags
  - gdcm::QueryBase, [614](#)
  - gdcm::QueryImage, [617](#)
  - gdcm::QueryPatient, [619](#)
  - gdcm::QuerySeries, [621](#)
  - gdcm::QueryStudy, [623](#)
- GetHighBit
  - gdcm::PixelFormat, [574](#)
- GetHostName
  - gdcm::System, [735](#)
- GetIE
  - gdcm::IODEntry, [461](#)
- GetIOD
  - gdcm::IODs, [463](#)
  - gdcm::SOPClassUIDToIOD, [690](#)
- GetIODEntry
  - gdcm::IOD, [460](#)
- GetIODFromFile

- gdcmm::Defs, 307
- GetIODFromSOPClassUID
  - gdcmm::SOPClassUIDToIOD, 690
- GetIODNameFromMediaStorage
  - gdcmm::Defs, 307
- GetIODs
  - gdcmm::Defs, 307, 308
- GetIconImage
  - gdcmm::IconImageFilter, 409
  - gdcmm::IconImageGenerator, 411
  - gdcmm::Pixmap, 579
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
- GetIconImagePort
  - vtkGDCMImageReader2, 868
- GetImage
  - gdcmm::ImageReader, 445
  - gdcmm::ImageWriter, 453
  - gdcmm::PixmapWriter, 587
  - gdcmm::SplitMosaicFilter, 697
- GetImplementationClassUID
  - gdcmm::FileMetaInformation, 382
- GetImplementationVersionName
  - gdcmm::FileMetaInformation, 382
- GetIndex
  - gdcmm::SwapCode, 731
  - gdcmm::VM, 849
- GetInitialized
  - gdcmm::CAPICryptographicMessageSyntax, 232
- GetInput
  - gdcmm::ImageToImageFilter, 451
  - gdcmm::PixmapToPixmapFilter, 585
  - vtkImageColorViewer, 895
- GetInputFilename
  - gdcmm::DictConverter, 317
- GetInstance
  - gdcmm::Global, 405
- GetIntercept
  - gdcmm::Image, 415
  - gdcmm::Rescaler, 636
- GetInterfile
  - gdcmm::CSAHeader, 272
- GetInternal
  - gdcmm::Preamble, 592
- GetIsCommand
  - gdcmm::network::PresentationDataValue, 600
- GetIsLastFragment
  - gdcmm::network::PresentationDataValue, 600
- GetItem
  - gdcmm::SequenceOfItems, 668
- GetKey
  - gdcmm::CSAElement, 266
- GetKeys
  - gdcmm::Scanner, 646
- GetKeyword
  - gdcmm::DictEntry, 319
- GetKeywordFromTag
  - gdcmm::Dict, 315
- GetLUT
  - gdcmm::Bitmap, 213
  - gdcmm::ImageCodec, 433
  - gdcmm::ImageHelper, 441
  - gdcmm::LookupTable, 499
- GetLUTDescriptor
  - gdcmm::LookupTable, 499
- GetLUTLength
  - gdcmm::LookupTable, 499
- GetLabel
  - gdcmm::Orientation, 545
- GetLastElement
  - gdcmm::ParseException, 553
- GetLastSystemError
  - gdcmm::System, 735
- GetLength
  - gdcmm::ByteValue, 227
  - gdcmm::CP246ExplicitDataElement, 259
  - gdcmm::DataElement, 286
  - gdcmm::DataSet, 299
  - gdcmm::Element, 337
  - gdcmm::Element< TVR, VM::VM1\_n >, 340
  - gdcmm::Element< VR::AS, VM::VM5 >, 348
  - gdcmm::ExplicitDataElement, 363
  - gdcmm::ExplicitImplicitDataElement, 365
  - gdcmm::Fragment, 403
  - gdcmm::ImplicitDataElement, 457
  - gdcmm::Item, 469
  - gdcmm::Preamble, 592
  - gdcmm::SequenceOfFragments, 662
  - gdcmm::SequenceOfItems, 668
  - gdcmm::Tag, 744
  - gdcmm::UNExplicitDataElement, 831
  - gdcmm::UNExplicitImplicitDataElement, 833
  - gdcmm::VL, 845
  - gdcmm::VM, 849
  - gdcmm::VR, 854
  - gdcmm::VR16ExplicitDataElement, 856
  - gdcmm::Value, 842
- GetLocaleCharset
  - gdcmm::System, 735
- GetLossless
  - gdcmm::JPEGCodec, 485
  - gdcmm::JPEGLSCodec, 489
- GetLossyFlag
  - gdcmm::ImageCodec, 433
- GetLossyFlagFromFile
  - gdcmm::Testing, 752
- GetMD5DataImage
  - gdcmm::Testing, 752



- GetMD5DataImages
  - gdcm::Testing, [752](#)
- GetMD5FromBrokenFile
  - gdcm::Testing, [752](#)
- GetMD5FromFile
  - gdcm::Testing, [752](#)
- GetMD5MetaImage
  - vtkGDCMTesting, [885](#)
- GetMHDMD5FromFile
  - vtkGDCMTesting, [885](#)
- GetMPType
  - gdcm::MeshPrimitive, [519](#)
- GetMPTypeString
  - gdcm::MeshPrimitive, [519](#)
- GetMRImageSeriesUIDs
  - gdcm::DirectoryHelper, [331](#)
- GetMSString
  - gdcm::MediaStorage, [511](#)
- GetMSType
  - gdcm::MediaStorage, [511](#)
- GetMTime
  - vtkImageMapToColors16, [901](#)
- GetMacro
  - gdcm::Macros, [503](#)
- GetMacroEntry
  - gdcm::Macro, [502](#)
- GetMacros
  - gdcm::Defs, [308](#)
- GetMajorAxisFromPatientRelativeDirectionCosine
  - gdcm::Orientation, [545](#)
- GetMajorVersion
  - gdcm::Version, [844](#)
- GetManifold
  - gdcm::Surface, [720](#)
- GetMapping
  - gdcm::Scanner, [646](#)
- GetMappingFromTagToValue
  - gdcm::Scanner, [646](#)
- GetMappings
  - gdcm::Scanner, [646](#)
- GetMax
  - gdcm::PixelFormat, [574](#)
- GetMaxLength
  - gdcm::PersonName, [566](#)
- GetMaxPDULength
  - gdcm::network::ULConnectionInfo, [823](#)
- GetMaxPDUSize
  - gdcm::network::ULConnection, [820](#)
- GetMaximumLength
  - gdcm::network::MaximumLengthSub, [504](#)
- GetMaximumLengthSub
  - gdcm::network::UserInformation, [837](#)
- GetMaximumPointDistance
  - gdcm::Surface, [720](#)
- GetMeanPointDistance
  - gdcm::Surface, [720](#)
- GetMediaStorage
  - gdcm::DataSet, [299](#)
  - gdcm::FileMetaInformation, [382](#)
- GetMediaStorageAsString
  - gdcm::FileMetaInformation, [382](#)
- GetMediaStorageDataFile
  - gdcm::Testing, [752](#)
- GetMediaStorageDataFiles
  - gdcm::Testing, [752](#)
- GetMediaStorageFromFile
  - gdcm::Testing, [752](#)
- GetMeshPrimitive
  - gdcm::Surface, [720](#)
- GetMessageHeader
  - gdcm::network::PresentationDataValue, [600](#)
- GetMetaInformationTS
  - gdcm::FileMetaInformation, [382](#)
- GetMin
  - gdcm::PixelFormat, [574](#)
- GetMinorVersion
  - gdcm::Version, [844](#)
- GetModality
  - gdcm::MediaStorage, [511](#)
- GetModalityDimension
  - gdcm::MediaStorage, [511](#)
- GetModule
  - gdcm::Modules, [526](#)
- GetModuleEntry
  - gdcm::NestedModuleEntries, [533](#)
- GetModuleEntryInMacros
  - gdcm::Module, [522](#)
- GetModules
  - gdcm::Defs, [308](#)
- GetName
  - gdcm::CSAElement, [266](#)
  - gdcm::CSAHeaderDictEntry, [275](#)
  - gdcm::DictEntry, [319](#)
  - gdcm::Filename, [385](#)
  - gdcm::GroupDict, [407](#)
  - gdcm::IODEntry, [461](#)
  - gdcm::Macro, [502](#)
  - gdcm::Module, [522](#)
  - gdcm::ModuleEntry, [524](#)
  - gdcm::PDBelement, [560](#)
  - gdcm::QueryBase, [614](#)
  - gdcm::QueryImage, [618](#)
  - gdcm::QueryPatient, [620](#)
  - gdcm::QuerySeries, [622](#)
  - gdcm::QueryStudy, [624](#)
  - gdcm::UIDs, [785](#)
  - gdcm::network::AbstractSyntax, [151](#)
  - gdcm::network::ApplicationContext, [161](#)

- gdcmm::network::TransferSyntaxSub, [762](#)
- GetNeedByteSwap
  - gdcmm::Bitmap, [213](#)
  - gdcmm::ImageCodec, [433](#)
- GetNegotiatedType
  - gdcmm::TransferSyntax, [760](#)
- GetNestedDataSet
  - gdcmm::Item, [469](#)
- GetNextSingleSerieUIDFileSet
  - gdcmm::SerieHelper, [672](#)
- GetNoOfItems
  - gdcmm::CSAElement, [266](#)
- GetNumberOfComponents
  - gdcmm::PersonName, [566](#)
- GetNumberOfContourReferencedFrameOfReferences
  - vtkRTStructSetProperties, [914](#)
- GetNumberOfCurves
  - gdcmm::Curve, [282](#)
  - gdcmm::Pixmap, [579](#)
- GetNumberOfDimensions
  - gdcmm::Bitmap, [213](#)
  - gdcmm::ImageCodec, [433](#)
- GetNumberOfElementsFromArray
  - gdcmm::VM, [849](#)
- GetNumberOfFileNames
  - gdcmm::Testing, [752](#)
- GetNumberOfFilenames
  - gdcmm::FilenameGenerator, [390](#)
- GetNumberOfFragments
  - gdcmm::SequenceOfFragments, [662](#)
- GetNumberOfIODs
  - gdcmm::IOD, [460](#)
- GetNumberOfIconImages
  - gdcmm::IconImageFilter, [410](#)
- GetNumberOfItems
  - gdcmm::SequenceOfItems, [668](#)
- GetNumberOfMD5DataImages
  - gdcmm::Testing, [752](#)
- GetNumberOfMD5MetalImages
  - vtkGDCMTesting, [885](#)
- GetNumberOfMSString
  - gdcmm::MediaStorage, [511](#)
- GetNumberOfMSType
  - gdcmm::MediaStorage, [511](#)
- GetNumberOfMediaStorageDataFiles
  - gdcmm::Testing, [752](#)
- GetNumberOfModality
  - gdcmm::MediaStorage, [511](#)
- GetNumberOfModuleEntries
  - gdcmm::NestedModuleEntries, [533](#)
- GetNumberOfOverlays
  - gdcmm::Pixmap, [579](#)
- GetNumberOfPoints
  - gdcmm::Curve, [282](#)
- GetNumberOfPresentationContext
  - gdcmm::network::AAssociateRQPDU, [148](#)
- GetNumberOfPresentationContextAC
  - gdcmm::network::AAssociateACPDU, [143](#)
- GetNumberOfPresentationDataValues
  - gdcmm::network::PDataTFPDU, [558](#)
- GetNumberOfPrimitivesData
  - gdcmm::MeshPrimitive, [519](#)
- GetNumberOfReferencedFrameOfReferences
  - vtkRTStructSetProperties, [914](#)
- GetNumberOfSOPClassToIOD
  - gdcmm::SOPClassUIDToIOD, [690](#)
- GetNumberOfSegments
  - gdcmm::SegmentWriter, [658](#)
- GetNumberOfStructureSetROIs
  - vtkRTStructSetProperties, [914](#)
- GetNumberOfSurfacePoints
  - gdcmm::Surface, [720](#)
- GetNumberOfSurfaces
  - gdcmm::SurfaceReader, [727](#)
  - gdcmm::SurfaceWriter, [729](#)
- GetNumberOfTransferSyntaxStrings
  - gdcmm::UIDs, [785](#)
- GetNumberOfTransferSyntaxes
  - gdcmm::PresentationContext, [593](#)
  - gdcmm::network::PresentationContextRQ, [598](#)
- GetNumberOfValues
  - gdcmm::Attribute, [172](#)
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
- GetNumberOfVectors
  - gdcmm::Surface, [720](#)
- GetObliquityThresholdCosineValue
  - gdcmm::Orientation, [546](#)
- GetOffScreenRendering
  - vtkImageColorViewer, [895](#)
- GetOptionalTags
  - gdcmm::QueryBase, [615](#)
  - gdcmm::QueryImage, [618](#)
  - gdcmm::QueryPatient, [620](#)
  - gdcmm::QuerySeries, [622](#)
  - gdcmm::QueryStudy, [624](#)
- GetOrderedValues
  - gdcmm::Scanner, [646](#)
- GetOrigin
  - gdcmm::Image, [415](#)
  - gdcmm::Overlay, [550](#)
- GetOriginValue
  - gdcmm::ImageHelper, [441](#)
- GetOutput
  - gdcmm::ImageConverter, [436](#)
- GetOutput

- gdcm::BitmapToBitmapFilter, 218
- gdcm::ImageToImageFilter, 451
- gdcm::PixmapToPixmapFilter, 585
- GetOutputAsBitmap
  - gdcm::BitmapToBitmapFilter, 218
- GetOutputAsPixmap
  - gdcm::PixmapToPixmapFilter, 585
- GetOutputFilename
  - gdcm::DictConverter, 317
- GetOutputType
  - gdcm::DictConverter, 317
- GetOverlay
  - gdcm::Pixmap, 579
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
- GetOverlayData
  - gdcm::Overlay, 550
- GetOverlayPort
  - vtkGDCMImageReader2, 868
- GetOverlayTypeAsString
  - gdcm::Overlay, 550
- GetOverlayTypeFromString
  - gdcm::Overlay, 550
- GetOverlayVisibility
  - vtkImageColorViewer, 896
- GetOwner
  - gdcm::PrivateTag, 607
- GetPDBEEnd
  - gdcm::PDBHeader, 562
- GetPDBElementByName
  - gdcm::PDBHeader, 562
- GetPDBInfoTag
  - gdcm::PDBHeader, 562
- GetPDUs
  - gdcm::network::ULEvent, 827
- GetPDVs
  - gdcm::network::PDUFactory, 565
- GetPIString
  - gdcm::PhotometricInterpretation, 571
- GetPIType
  - gdcm::PhotometricInterpretation, 571
- GetPath
  - gdcm::Filename, 385
- GetPattern
  - gdcm::FilenameGenerator, 390
- GetPermissions
  - gdcm::System, 735
- GetPhotometricInterpretation
  - gdcm::Bitmap, 213
  - gdcm::ImageChangePhotometricInterpretation, 421
  - gdcm::ImageCodec, 433
- GetPhotometricInterpretationValue
  - gdcm::ImageHelper, 441
- GetPixelFormat
  - gdcm::Bitmap, 213, 214
  - gdcm::ImageCodec, 433
- GetPixelFormatValue
  - gdcm::ImageHelper, 441
- GetPixelRepresentation
  - gdcm::PixelFormat, 575
- GetPixelSize
  - gdcm::PixelFormat, 575
- GetPixelSpacingDataRoot
  - gdcm::Testing, 752
- GetPixmap
  - gdcm::IconImageGenerator, 412
  - gdcm::PixmapReader, 582
  - gdcm::PixmapWriter, 587
- GetPlanarConfiguration
  - gdcm::Bitmap, 214
  - gdcm::ImageChangePlanarConfiguration, 424
  - gdcm::ImageCodec, 433
- GetPlanarConfigurationValue
  - gdcm::ImageHelper, 441
- GetPointCoordinatesData
  - gdcm::Surface, 720
- GetPointPositionAccuracy
  - gdcm::Surface, 720
- GetPointer
  - gdcm::ByteValue, 227
  - gdcm::LookupTable, 499
  - gdcm::SmartPointer, 687
  - vtkLookupTable16, 911
- GetPointerFromElement
  - gdcm::ImageHelper, 441
- GetPointsBoundingBoxCoordinates
  - gdcm::Surface, 720
- GetPosition
  - vtkImageColorViewer, 896
- GetPreamble
  - gdcm::FileMetaInformation, 382
- GetPrefix
  - gdcm::FilenameGenerator, 390
- GetPresentationContext
  - gdcm::network::AAssociateRQPDU, 148
- GetPresentationContextAC
  - gdcm::network::AAssociateACPDU, 143
- GetPresentationContextACByID
  - gdcm::network::ULConnection, 820
- GetPresentationContextByAbstractSyntax
  - gdcm::network::AAssociateRQPDU, 148
- GetPresentationContextByID
  - gdcm::network::AAssociateRQPDU, 148
- GetPresentationContextID
  - gdcm::PresentationContext, 593
  - gdcm::network::PresentationContextAC, 594
  - gdcm::network::PresentationContextRQ, 598
  - gdcm::network::PresentationDataValue, 600

- GetPresentationContextIDFromPresentationContext
  - gdcm::network::ULConnection, [820](#)
- GetPresentationContextRQByID
  - gdcm::network::ULConnection, [820](#)
- GetPresentationContexts
  - gdcm::PresentationContextGenerator, [597](#)
  - gdcm::network::AAssociateRQPDU, [148](#)
  - gdcm::network::ULConnection, [820](#)
- GetPresentationDataValue
  - gdcm::network::PDataTFPDU, [558](#)
- GetPrettyPrint
  - gdcm::JSON, [491](#)
- GetPrimitiveData
  - gdcm::MeshPrimitive, [519](#)
- GetPrimitiveType
  - gdcm::MeshPrimitive, [519](#)
- GetPrimitivesData
  - gdcm::MeshPrimitive, [519](#)
- GetPrintStyle
  - gdcm::Printer, [603](#)
  - gdcm::XMLPrinter, [925](#)
- GetPrivateCreator
  - gdcm::DataSet, [299](#)
  - gdcm::Tag, [744](#)
- GetPrivateDict
  - gdcm::Dicts, [324](#)
  - gdcm::XMLPrivateDictReader, [927](#)
- GetProcessingAlgorithm
  - gdcm::Surface, [720](#)
- GetProgress
  - gdcm::ProgressEvent, [609](#)
- GetPropertyCategory
  - gdcm::Segment, [650](#)
- GetPropertyType
  - gdcm::Segment, [650](#)
- GetProtocol
  - gdcm::network::ULConnection, [820](#)
- GetPublicDict
  - gdcm::Dicts, [324](#)
- GetQuality
  - gdcm::JPEG2000Codec, [479](#)
  - gdcm::JPEGCodec, [485](#)
- GetQueryDataSet
  - gdcm::BaseRootQuery, [203](#)
- GetQueryLevel
  - gdcm::QueryBase, [615](#)
  - gdcm::QueryImage, [618](#)
  - gdcm::QueryPatient, [620](#)
  - gdcm::QuerySeries, [622](#)
  - gdcm::QueryStudy, [624](#)
- GetQueryLevelFromQueryRoot
  - gdcm::BaseRootQuery, [203](#)
- GetQueryLevelFromString
  - gdcm::BaseRootQuery, [203](#)
- GetQueryLevelString
  - gdcm::BaseRootQuery, [203](#)
- GetRAWMD5FromFile
  - vtkGDCMTesting, [885](#)
- GetRTStructSeriesUIDs
  - gdcm::DirectoryHelper, [331](#)
- GetRate
  - gdcm::JPEG2000Codec, [479](#)
- GetReason
  - gdcm::network::PresentationContextAC, [595](#)
- GetRecommendedDisplayCIELabValue
  - gdcm::Surface, [720](#)
- GetRecommendedDisplayGrayscaleValue
  - gdcm::Surface, [720](#)
- GetRecommendedPresentationOpacity
  - gdcm::Surface, [720](#)
- GetRecommendedPresentationType
  - gdcm::Surface, [720](#)
- GetRef
  - gdcm::IODEntry, [461](#)
- GetReferencedFrameOfReferenceClassUID
  - vtkRTStructSetProperties, [914](#)
- GetReferencedFrameOfReferenceInstanceUID
  - vtkRTStructSetProperties, [914](#)
- GetRegion
  - gdcm::ImageRegionReader, [448](#)
- GetRequiredTags
  - gdcm::QueryBase, [615](#)
  - gdcm::QueryImage, [618](#)
  - gdcm::QueryPatient, [620](#)
  - gdcm::QuerySeries, [622](#)
  - gdcm::QueryStudy, [624](#)
- GetRescaleInterceptSlopeValue
  - gdcm::ImageHelper, [441](#)
- GetReserved43\_74
  - gdcm::network::AAssociateRQPDU, [148](#)
- GetResponses
  - gdcm::network::ULBasicCallback, [818](#)
- GetRetired
  - gdcm::DictEntry, [319](#)
- GetRoot
  - gdcm::UIDGenerator, [767](#)
- GetRows
  - gdcm::Bitmap, [214](#)
  - gdcm::Overlay, [550](#)
- GetSOPClassUID
  - gdcm::DirectoryHelper, [332](#)
- GetSOPClassUIDFromIOD
  - gdcm::SOPClassUIDToIOD, [690](#)
- GetSOPClassUIDToIOD
  - gdcm::SOPClassUIDToIOD, [690](#)
- GetSOPClassUIDToIODs
  - gdcm::SOPClassUIDToIOD, [690](#)
- GetSTATES

- gdcm::Surface, [721](#)
- GetSTATESString
  - gdcm::Surface, [721](#)
- GetSamplesPerPixel
  - gdcm::PhotometricInterpretation, [571](#)
  - gdcm::PixelFormat, [575](#)
- GetScalarType
  - gdcm::PixelFormat, [575](#)
- GetScalarTypeAsString
  - gdcm::PixelFormat, [575](#)
- GetScanner
  - gdcm::DICOMDIRGenerator, [313](#)
- GetSegment
  - gdcm::SegmentWriter, [658](#)
- GetSegmentAlgorithmName
  - gdcm::Segment, [650](#)
- GetSegmentAlgorithmType
  - gdcm::Segment, [651](#)
- GetSegmentDescription
  - gdcm::Segment, [651](#)
- GetSegmentLabel
  - gdcm::Segment, [651](#)
- GetSegmentNumber
  - gdcm::Segment, [651](#)
- GetSegments
  - gdcm::SegmentReader, [656](#)
  - gdcm::SegmentWriter, [658](#)
- GetSelectedPrivateGroupOffsetFromFile
  - gdcm::Testing, [752](#)
- GetSelectedTagsOffsetFromFile
  - gdcm::Testing, [753](#)
- GetSequenceOfFragments
  - gdcm::DataElement, [286](#)
- GetSeriesUIDsBySOPClassUID
  - gdcm::DirectoryHelper, [332](#)
- GetSize
  - gdcm::VR, [854](#)
  - vtkImageColorViewer, [896](#)
- GetSizeof
  - gdcm::VR, [854](#)
- GetSliceMax
  - vtkImageColorViewer, [896](#)
- GetSliceMin
  - vtkImageColorViewer, [896](#)
- GetSliceRange
  - vtkImageColorViewer, [896](#)
- GetSlope
  - gdcm::Image, [416](#)
  - gdcm::Rescaler, [636](#)
- GetSourceApplicationEntityTitle
  - gdcm::FileMetaInformation, [382](#)
- GetSourceDirectory
  - gdcm::Testing, [753](#)
- GetSpacing
  - gdcm::Image, [416](#)
- GetSpacingTagFromMediaStorage
  - gdcm::ImageHelper, [442](#)
- GetSpacingValue
  - gdcm::ImageHelper, [442](#)
- GetStart
  - gdcm::ByteBuffer, [222](#)
- GetState
  - gdcm::network::ULConnection, [820](#)
- GetStateIndex
  - gdcm::network, [136](#)
- GetStream
  - gdcm::Trace, [756](#)
- GetStreamCurrentPosition
  - gdcm::Reader, [630](#)
- GetStreamOffsetFromFile
  - gdcm::Testing, [753](#)
- GetStreamPtr
  - gdcm::Reader, [630](#)
  - gdcm::Writer, [920](#)
- GetString
  - gdcm::MediaStorage, [511](#)
  - gdcm::PhotometricInterpretation, [571](#)
  - gdcm::TransferSyntax, [760](#)
  - gdcm::UIDs, [786](#)
- GetStringValueFromTag
  - gdcm::DirectoryHelper, [332](#)
- GetStructureSetObservationNumber
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROIDescription
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROIGenerationAlgorithm
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROIName
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROINumber
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROIObservationLabel
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetROIRefFrameRefUID
  - vtkRTStructSetProperties, [914](#)
- GetStructureSetRTROIInterpretedType
  - vtkRTStructSetProperties, [914](#)
- GetSurface
  - gdcm::Segment, [651](#)
- GetSurfaceComments
  - gdcm::Surface, [721](#)
- GetSurfaceCount
  - gdcm::Segment, [651](#)
- GetSurfaceNumber
  - gdcm::Surface, [721](#)
- GetSurfaceProcessing
  - gdcm::Surface, [721](#)
- GetSurfaceProcessingDescription

- gdcm::Surface, [721](#)
- GetSurfaceProcessingRatio
  - gdcm::Surface, [721](#)
- GetSurfaces
  - gdcm::Segment, [651](#)
- GetSwapCode
  - gdcm::TransferSyntax, [760](#)
- GetSwapCodeString
  - gdcm::SwapCode, [731](#)
- GetSyngoDT
  - gdcm::CSAElement, [266](#)
- GetTSString
  - gdcm::TransferSyntax, [760](#)
- GetTSType
  - gdcm::TransferSyntax, [760](#)
- GetTable
  - gdcm::SequenceOfFragments, [662](#), [663](#)
- GetTableEntry
  - gdcm::Table, [737](#)
- GetTag
  - gdcm::AnonymizeEvent, [153](#)
  - gdcm::Attribute, [172](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::DataElement, [287](#)
- GetTagListByLevel
  - gdcm::BaseRootQuery, [203](#)
  - gdcm::FindPatientRootQuery, [398](#)
  - gdcm::FindStudyRootQuery, [401](#)
  - gdcm::MovePatientRootQuery, [528](#)
  - gdcm::MoveStudyRootQuery, [530](#)
- GetTempDirectory
  - gdcm::Testing, [753](#)
- GetTempDirectoryW
  - gdcm::Testing, [753](#)
- GetTempFilename
  - gdcm::Testing, [753](#)
- GetTempFilenameW
  - gdcm::Testing, [753](#)
- GetTimeout
  - gdcm::ServiceClassUser, [677](#)
  - gdcm::network::ARTIMTimer, [167](#)
- GetTimer
  - gdcm::network::ULConnection, [820](#)
- GetTimezoneOffsetFromUTC
  - gdcm::System, [735](#)
- GetToplevel
  - gdcm::Directory, [330](#)
- GetTransferSyntax
  - gdcm::Bitmap, [214](#)
  - gdcm::ImageChangeTransferSyntax, [427](#)
  - gdcm::PresentationContext, [594](#)
  - gdcm::network::PresentationContextAC, [595](#)
  - gdcm::network::PresentationContextRQ, [598](#), [599](#)
- GetTransferSyntaxString
  - gdcm::UIDs, [786](#)
- GetTransferSyntaxStrings
  - gdcm::UIDs, [786](#)
- GetTransferSyntaxes
  - gdcm::network::PresentationContextRQ, [599](#)
- GetType
  - gdcm::ModuleEntry, [524](#)
  - gdcm::Orientation, [546](#)
  - gdcm::Overlay, [550](#)
  - gdcm::PhotometricInterpretation, [571](#)
- GetTypeAsEnum
  - gdcm::Overlay, [550](#)
- GetTypeFromTag
  - gdcm::Defs, [308](#)
  - gdcm::IOD, [460](#)
- GetTypeOfData
  - gdcm::Curve, [282](#)
- GetTypeOfDataDescription
  - gdcm::Curve, [282](#)
- GetTypeString
  - gdcm::Type, [765](#)
- GetTypeType
  - gdcm::Type, [765](#)
- GetUIDName
  - gdcm::UIDs, [786](#)
- GetUIDString
  - gdcm::UIDs, [786](#)
- GetUniqueTags
  - gdcm::QueryBase, [615](#)
  - gdcm::QueryImage, [618](#)
  - gdcm::QueryPatient, [620](#)
  - gdcm::QuerySeries, [622](#)
  - gdcm::QueryStudy, [624](#)
- GetUnpackBuffer
  - gdcm::Overlay, [550](#)
- GetUnpackBufferLength
  - gdcm::Overlay, [550](#)
- GetUsage
  - gdcm::IODEntry, [461](#)
- GetUsageString
  - gdcm::Usage, [835](#)
- GetUsageType
  - gdcm::IODEntry, [462](#)
  - gdcm::Usage, [835](#)
- GetUserData
  - gdcm::Parser, [555](#)
- GetUserInformation
  - gdcm::network::AAssociateACPDU, [143](#)
  - gdcm::network::AAssociateRQPDU, [148](#)
- GetVIEWType
  - gdcm::Surface, [721](#)



- GetVIEWTypeString
  - gdcm::Surface, [721](#)
- GetVL
  - gdcm::DataElement, [287](#)
- GetVL16Max
  - gdcm::VL, [845](#)
- GetVL32Max
  - gdcm::VL, [845](#)
- GetVM
  - gdcm::Attribute, [173](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_3 >, [181](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_8 >, [182](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM2↔\_2n >, [188](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM2\_n >, [189](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM3↔\_3n >, [191](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM3\_n >, [192](#)
  - gdcm::CSAElement, [267](#)
  - gdcm::CSAHeaderDictEntry, [275](#)
  - gdcm::DictEntry, [319](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [340](#)
- GetVMString
  - gdcm::VM, [849](#)
- GetVMType
  - gdcm::VM, [849](#)
- GetVMTypeFromLength
  - gdcm::VM, [850](#)
- GetVR
  - gdcm::Attribute, [173](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::CSAElement, [267](#)
  - gdcm::CSAHeaderDictEntry, [275](#)
  - gdcm::DataElement, [288](#)
  - gdcm::DictEntry, [319](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [340](#)
- GetVRFromTag
  - gdcm, [126](#)
- GetVRString
  - gdcm::VR, [854](#)
- GetVRStringFromFile
  - gdcm::VR, [854](#)
- GetVRType
  - gdcm::VR, [854](#)
- GetVRTypeFromFile
  - gdcm::VR, [854](#)
- GetVTKDataRoot
  - vtkGDCMTesting, [885](#)
- GetValidatedFile
  - gdcm::Validate, [840](#)
- GetValue
  - gdcm::Attribute, [172](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::CSAElement, [266](#)
  - gdcm::DataElement, [287](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [340](#)
  - gdcm::PDSElement, [560](#)
  - gdcm::Scanner, [646](#)
- GetValueAsSQ
  - gdcm::DataElement, [287](#)
- GetValues
  - gdcm::Attribute, [172](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [178](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::Element, [337](#)
  - gdcm::Scanner, [647](#)
- GetVectorAccuracy
  - gdcm::Surface, [721](#)
- GetVectorCoordinateData
  - gdcm::Surface, [721](#)
- GetVectorDimensionality
  - gdcm::Surface, [721](#)
- GetVersion
  - gdcm::Version, [844](#)
- GetWarningFlag
  - gdcm::Trace, [756](#)
- GetWarningStream
  - gdcm::Trace, [756](#)
- GetWindowName
  - vtkImageColorViewer, [896](#)
- GetXMax
  - gdcm::BoxRegion, [221](#)
- GetXMin
  - gdcm::BoxRegion, [221](#)
- GetYMax
  - gdcm::BoxRegion, [221](#)
- GetYMin
  - gdcm::BoxRegion, [221](#)
- GetZMax

- gdcmm::BoxRegion, 221
- GetZMin
  - gdcmm::BoxRegion, 221
- GetZSpacing
  - gdcmm::IPPSorter, 465
- GetZSpacingTagFromMediaStorage
  - gdcmm::ImageHelper, 442
- GetZSpacingTolerance
  - gdcmm::IPPSorter, 466
- Global
  - gdcmm::Defs, 308
  - gdcmm::Dicts, 324
  - gdcmm::Global, 405
- GlobalInstance
  - gdcmm, 131
- GrabOverlayFromPixelData
  - gdcmm::Overlay, 550
- Graphics
  - gdcmm::Overlay, 549
- GrayscaleSoftcopyPresentationStateStorageSOPClass
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 776
- green
  - gdcmm::terminal, 137
- group
  - gdcmm::SerieHelper::Rule, 642
- GroupDict
  - gdcmm::GroupDict, 407
- GroupStringVector
  - gdcmm::GroupDict, 407
- GuessFromModality
  - gdcmm::MediaStorage, 511
- HSV
  - gdcmm::PhotometricInterpretation, 570
- HandleBulkData
  - gdcmm::XMLPrinter, 925
- HandleDataSet
  - gdcmm::network::ULBasicCallback, 818
  - gdcmm::network::ULConnectionCallback, 822
  - gdcmm::network::ULWritingCallback, 829
- HandleDescription
  - gdcmm::XMLDictReader, 923
  - gdcmm::XMLPrivateDictReader, 927
- HandleEntry
  - gdcmm::XMLDictReader, 923
  - gdcmm::XMLPrivateDictReader, 927
- HandleEvent
  - gdcmm::network::ULTransitionTable, 828
- HandleIOD
  - gdcmm::TableReader, 739
- HandleIODEntry
  - gdcmm::TableReader, 739
- HandleMacro
  - gdcmm::TableReader, 739
- HandleMacroEntry
  - gdcmm::TableReader, 739
- HandleMacroEntryDescription
  - gdcmm::TableReader, 739
- HandleModule
  - gdcmm::TableReader, 740
- HandleModuleEntry
  - gdcmm::TableReader, 740
- HandleModuleEntryDescription
  - gdcmm::TableReader, 740
- HandleModuleInclude
  - gdcmm::TableReader, 740
- HandleResponse
  - gdcmm::network::ULBasicCallback, 818
  - gdcmm::network::ULConnectionCallback, 822
  - gdcmm::network::ULWritingCallback, 829
- HangingProtocolInformationModelFIND
  - gdcmm::UIDs, 778
- HangingProtocolInformationModelMOVE
  - gdcmm::UIDs, 778
- HangingProtocolStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 778
- HardcopyColorImageStorageSOPClassRetired
  - gdcmm::UIDs, 775
- HardcopyGrayscaleImageStorage
  - gdcmm::MediaStorage, 509
- HardcopyGrayscaleImageStorageSOPClassRetired
  - gdcmm::UIDs, 775
- HasObserver
  - gdcmm::Subject, 715
- HemodynamicWaveformStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 775
- hidden
  - gdcmm::terminal, 137
- ICBM452T1FrameofReference
  - gdcmm::UIDs, 774
- ICBMSingleSubjectMRIFrameofReference
  - gdcmm::UIDs, 774
- INT12
  - gdcmm::PixelFormat, 574
- INT16
  - gdcmm::PixelFormat, 574
- INT32
  - gdcmm::PixelFormat, 574
- INT64
  - gdcmm::PixelFormat, 574
- INT8
  - gdcmm::PixelFormat, 573
- INTERFILE
  - gdcmm::CSAHeader, 270



- INVALID
  - gdcm::VR, [852](#)
- IOD
  - gdcm::IOD, [459](#)
- IODEntry
  - gdcm::IODEntry, [461](#)
- IODMapType
  - gdcm::IODs, [463](#)
- IODMapTypeConstIterator
  - gdcm::IODs, [463](#)
- IODName
  - gdcm::IODs, [463](#)
- IODs
  - gdcm::IODs, [463](#)
- IPPSorter
  - gdcm::IPPSorter, [465](#)
- IS
  - gdcm::VR, [853](#)
- Icon
  - gdcm::Pixmap, [579](#)
- IconDataScalarType
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- IconImage
  - gdcm, [124](#)
- IconImageDataExtent
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- IconImageFilter
  - gdcm::IconImageFilter, [409](#)
- IconImageGenerator
  - gdcm::IconImageGenerator, [411](#)
- IconNumberOfScalarComponents
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- ignore\_char
  - gdcm::ignore\_char, [413](#)
- Image
  - gdcm::Image, [415](#)
- ImageActor
  - vtkImageColorViewer, [898](#)
- ImageApplyLookupTable
  - gdcm::ImageApplyLookupTable, [419](#)
- ImageChangePhotometricInterpretation
  - gdcm::ImageChangePhotometricInterpretation, [421](#)
  - gdcm::ImageCodec, [434](#)
- ImageChangePlanarConfiguration
  - gdcm::ImageChangePlanarConfiguration, [424](#)
- ImageChangeTransferSyntax
  - gdcm::Bitmap, [216](#)
  - gdcm::ImageChangeTransferSyntax, [427](#)
- ImageCodec
  - gdcm::ImageCodec, [431](#)
- ImageConverter
  - gdcm::ImageConverter, [436](#)
- ImageFormat
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- ImageFragmentSplitter
  - gdcm::ImageFragmentSplitter, [439](#)
- ImageOrientationPatient
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- ImageOverlayBoxSOPClassRetired
  - gdcm::UIDs, [775](#)
- ImagePositionPatient
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
- ImagePositionPatientOrdering
  - gdcm::SerieHelper, [672](#)
- ImageReader
  - gdcm::ImageReader, [445](#)
- ImageRegionReader
  - gdcm::ImageRegionReader, [448](#)
  - gdcm::JPEG2000Codec, [479](#)
  - gdcm::JPEGCodec, [486](#)
  - gdcm::JPEGLSCodec, [490](#)
  - gdcm::RLECodec, [640](#)
- ImageToImageFilter
  - gdcm::ImageToImageFilter, [450](#)
- ImageWriter
  - gdcm::ImageWriter, [453](#)
- ImplementationClassUIDSub
  - gdcm::network::ImplementationClassUIDSub, [454](#)
- ImplementationUIDSub
  - gdcm::network::ImplementationUIDSub, [455](#)
- ImplementationVersionNameSub
  - gdcm::network::ImplementationVersionNameSub, [455](#)
- Implicit
  - gdcm::TransferSyntax, [759](#)
- ImplicitVRBigEndianACRNEMA
  - gdcm::TransferSyntax, [760](#)
- ImplicitVRBigEndianPrivateGE
  - gdcm::TransferSyntax, [759](#)
- ImplicitVRLittleEndian
  - gdcm::TransferSyntax, [759](#)
- ImplicitVRLittleEndianDefaultTransferSyntaxforDICOM
  - gdcm::UIDs, [772](#)
- IncompleteLUT
  - gdcm::LookupTable, [500](#)
- InitFromRQ
  - gdcm::network::AAssociateACPDU, [143](#)
- InitOpenSSL
  - gdcm::OpenSSLCryptoFactory, [538](#)
- Initialize
  - gdcm::network::ULConnectionInfo, [823](#)
- InitializeBlueLUT

- gdcmm::LookupTable, 499
- InitializeConnection
  - gdcmm::ServiceClassUser, 677
  - gdcmm::network::ULConnection, 820
- InitializeDataSet
  - gdcmm::BaseRootQuery, 203
  - gdcmm::FindPatientRootQuery, 399
  - gdcmm::FindStudyRootQuery, 401
  - gdcmm::MovePatientRootQuery, 528
  - gdcmm::MoveStudyRootQuery, 530
- InitializeGreenLUT
  - gdcmm::LookupTable, 499
- InitializeIncomingConnection
  - gdcmm::network::ULConnection, 820
- InitializeLUT
  - gdcmm::LookupTable, 499
- InitializeRTStructSet
  - vtkGDCMPolyDataWriter, 882
- InitializeRedLUT
  - gdcmm::LookupTable, 499
- Initialized
  - gdcmm::LookupTable, 499
- Input
  - gdcmm::BitmapToBitmapFilter, 218
- Insert
  - gdcmm::CommandDataSet, 252
  - gdcmm::DataSet, 299
  - gdcmm::FileMetaInformation, 382
  - gdcmm::GroupDict, 407
- InsertDataElement
  - gdcmm::DataSet, 299
  - gdcmm::Item, 470
- InsertEntry
  - gdcmm::Table, 737
- InstallPipeline
  - vtkImageColorViewer, 896
- InstanceAvailabilityNotificationSOPClass
  - gdcmm::UIDs, 777
- Interactor
  - vtkImageColorViewer, 898
- InteractorStyle
  - vtkImageColorViewer, 898
- Internal
  - gdcmm::ApplicationEntity, 162
  - gdcmm::Attribute, 175
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 180
  - gdcmm::Element, 337
  - gdcmm::Element< VR::AS, VM::VM5 >, 348
  - gdcmm::LookupTable, 500
  - gdcmm::UI, 766
- InternalCode
  - gdcmm::Coder, 245
  - gdcmm::JPEG12Codec, 473
  - gdcmm::JPEG16Codec, 475
  - gdcmm::JPEG8Codec, 481
- Internals
  - vtkRTStructSetProperties, 915
- Invalid
  - gdcmm::Overlay, 549
  - gdcmm::Usage, 835
- InverseRescale
  - gdcmm::Rescaler, 636
- InverseRescaleFunctionIntoBestFit
  - gdcmm::Rescaler, 636
- InvokeEvent
  - gdcmm::Subject, 715
- IsAETitleValid
  - gdcmm::network::AAssociateRQPDU, 148
- IsASCII
  - gdcmm::VR, 854
- IsASCII2
  - gdcmm::VR, 854
- IsBinary
  - gdcmm::VR, 854
- IsBinary2
  - gdcmm::VR, 854
- IsCompatible
  - gdcmm::PixelFormat, 575
- IsDual
  - gdcmm::VR, 854
- IsEmpty
  - gdcmm::Bitmap, 214
  - gdcmm::ByteValue, 228
  - gdcmm::CSAElement, 267
  - gdcmm::CSAHeaderDict, 274
  - gdcmm::Curve, 282
  - gdcmm::DataElement, 288
  - gdcmm::DataSet, 299
  - gdcmm::Defs, 308
  - gdcmm::Dict, 315
  - gdcmm::Dicts, 324
  - gdcmm::Filename, 385
  - gdcmm::Macros, 503
  - gdcmm::Modules, 526
  - gdcmm::Overlay, 550
  - gdcmm::Preamble, 592
  - gdcmm::PrivateDict, 605
  - gdcmm::SegmentHelper::BasicCodedEntry, 206
- IsEncapsulated
  - gdcmm::TransferSyntax, 760
- IsEncoded
  - gdcmm::TransferSyntax, 760
- IsExplicit
  - gdcmm::TransferSyntax, 761
- IsFrameEncoder
  - gdcmm::ImageCodec, 433
  - gdcmm::JPEG2000Codec, 479

- gdcm::JPEGCodec, 485
- gdcm::JPEGLSCodec, 489
- gdcm::RLECodec, 639
- IsGroupLength
  - gdcm::Tag, 744
- IsGroupXX
  - gdcm::Tag, 744
- IsIdentical
  - gdcm::Filename, 385
- IsIllegal
  - gdcm::Tag, 745
- IsImage
  - gdcm::MediaStorage, 511
- IsImplicit
  - gdcm::TransferSyntax, 761
- IsInPixelData
  - gdcm::Overlay, 550
- IsKey
  - gdcm::Scanner, 647
- IsLastFragment
  - gdcm::network::AAAbortPDU, 140
  - gdcm::network::AAAssociateACPDU, 143
  - gdcm::network::AAAssociateRJPDU, 145
  - gdcm::network::AAAssociateRQPDU, 148
  - gdcm::network::AReleaseRPPDU, 164
  - gdcm::network::AReleaseRQPDU, 166
  - gdcm::network::BasePDU, 200
  - gdcm::network::PDataTFPDU, 558
- IsLossless
  - gdcm::PhotometricInterpretation, 571
  - gdcm::TransferSyntax, 761
- IsLossy
  - gdcm::Bitmap, 214
  - gdcm::ImageCodec, 433
  - gdcm::PhotometricInterpretation, 571
  - gdcm::TransferSyntax, 761
- IsOdd
  - gdcm::VL, 845
- IsPresentationContextAccepted
  - gdcm::ServiceClassUser, 677
- IsPrintable
  - gdcm::ByteValue, 228
- IsPrivate
  - gdcm::Tag, 745
- IsPrivateCreator
  - gdcm::Tag, 745
- IsPublic
  - gdcm::Tag, 745
- IsRetired
  - gdcm::PhotometricInterpretation, 571
- IsRowEncoder
  - gdcm::ImageCodec, 433
  - gdcm::JPEG2000Codec, 479
  - gdcm::JPEGCodec, 485
  - gdcm::JPEGLSCodec, 489
  - gdcm::RLECodec, 640
- IsSameColorSpace
  - gdcm::PhotometricInterpretation, 571
- IsStateSuspension
  - gdcm::JPEG12Codec, 473
  - gdcm::JPEG16Codec, 475
  - gdcm::JPEG8Codec, 481
  - gdcm::JPEGCodec, 485
- IsSwap
  - gdcm::VR, 854
- IsTransferSyntaxCompatible
  - gdcm::Bitmap, 214
- IsUndefined
  - gdcm::MediaStorage, 512
  - gdcm::VL, 845
- IsUndefinedLength
  - gdcm::DataElement, 288
  - gdcm::SequenceOfItems, 668
- IsUnique
  - gdcm::DictEntry, 320
- IsVRFile
  - gdcm::VR, 854
- IsValid
  - gdcm::ApplicationEntity, 162
  - gdcm::BoxRegion, 221
  - gdcm::CodeString, 248
  - gdcm::DirectionCosines, 327
  - gdcm::FileMetaInformation, 382
  - gdcm::ImageCodec, 433
  - gdcm::JPEGCodec, 486
  - gdcm::LO, 495
  - gdcm::PixelFormat, 575
  - gdcm::Preamble, 592
  - gdcm::Region, 633
  - gdcm::String, 710
  - gdcm::TagPath, 749
  - gdcm::TransferSyntax, 761
  - gdcm::UIDGenerator, 767
  - gdcm::UUIDGenerator, 838
  - gdcm::VM, 850
  - gdcm::VR, 854
- IsZero
  - gdcm::Overlay, 551
- ItFileSetHt
  - gdcm::SerieHelper, 672
- Item
  - gdcm::Item, 469
- ItemVector
  - gdcm::SequenceOfItems, 667
- Items
  - gdcm::SequenceOfItems, 669
- Iterator
  - gdcm::CSAHeaderDict, 273

- gdcm::DataSet, [297](#)
- gdcm::Dict, [314](#)
- gdcm::SequenceOfFragments, [661](#)
- gdcm::SequenceOfItems, [667](#)
- iterator
  - gdcm::CodeString, [247](#)
  - gdcm::LO, [495](#)
  - gdcm::String, [709](#)
- JPEG12Codec
  - gdcm::JPEG12Codec, [473](#)
- JPEG16Codec
  - gdcm::JPEG16Codec, [475](#)
- JPEG2000
  - gdcm::TransferSyntax, [760](#)
- JPEG2000\_COMPRESSION
  - vtkGDCMImageWriter, [873](#)
- JPEG2000Codec
  - gdcm::JPEG2000Codec, [477](#)
- JPEG2000ImageCompression
  - gdcm::UIDs, [773](#)
- JPEG2000ImageCompressionLosslessOnly
  - gdcm::UIDs, [773](#)
- JPEG2000Lossless
  - gdcm::TransferSyntax, [760](#)
- JPEG2000Part2
  - gdcm::TransferSyntax, [760](#)
- JPEG2000Part2Lossless
  - gdcm::TransferSyntax, [760](#)
- JPEG2000Part2MulticomponentImageCompression
  - gdcm::UIDs, [773](#)
- JPEG2000Part2MulticomponentImageCompression↔
  - LosslessOnly
    - gdcm::UIDs, [773](#)
- JPEG8Codec
  - gdcm::JPEG8Codec, [481](#)
- JPEG\_COMPRESSION
  - vtkGDCMImageWriter, [873](#)
- JPEGBaselineProcess1
  - gdcm::TransferSyntax, [759](#)
- JPEGBaselineProcess1DefaultTransferSyntaxforLossyJ↔
  - PEG8BitImageCompression
    - gdcm::UIDs, [772](#)
- JPEGCodec
  - gdcm::JPEGCodec, [484](#)
- JPEGExtendedHierarchicalProcess1618Retired
  - gdcm::UIDs, [773](#)
- JPEGExtendedHierarchicalProcess1719Retired
  - gdcm::UIDs, [773](#)
- JPEGExtendedProcess24DefaultTransferSyntaxfor↔
  - LossyJPEG12BitImageCompressionProcess4only
    - gdcm::UIDs, [772](#)
- JPEGExtendedProcess2\_4
  - gdcm::TransferSyntax, [759](#)
- JPEGExtendedProcess35Retired
  - gdcm::UIDs, [772](#)
- JPEGExtendedProcess3\_5
  - gdcm::TransferSyntax, [759](#)
- JPEGFullProgressionHierarchicalProcess2426Retired
  - gdcm::UIDs, [773](#)
- JPEGFullProgressionHierarchicalProcess2527Retired
  - gdcm::UIDs, [773](#)
- JPEGFullProgressionNonHierarchicalProcess1012↔
  - Retired
    - gdcm::UIDs, [772](#)
- JPEGFullProgressionNonHierarchicalProcess1113↔
  - Retired
    - gdcm::UIDs, [772](#)
- JPEGFullProgressionProcess10\_12
  - gdcm::TransferSyntax, [759](#)
- JPEGLS\_COMPRESSION
  - vtkGDCMImageWriter, [873](#)
- JPEGLSCodec
  - gdcm::JPEGLSCodec, [488](#)
- JPEGLSLossless
  - gdcm::TransferSyntax, [759](#)
- JPEGLSLosslessImageCompression
  - gdcm::UIDs, [773](#)
- JPEGLSLossyNearLosslessImageCompression
  - gdcm::UIDs, [773](#)
- JPEGLSNearLossless
  - gdcm::TransferSyntax, [759](#)
- JPEGLosslessHierarchicalProcess28Retired
  - gdcm::UIDs, [773](#)
- JPEGLosslessHierarchicalProcess29Retired
  - gdcm::UIDs, [773](#)
- JPEGLosslessNonHierarchicalFirstOrderPrediction↔
  - Process14SelectionValue1DefaultTransfer↔
    - SyntaxforLosslessJPEGImageCompression
      - gdcm::UIDs, [773](#)
- JPEGLosslessNonHierarchicalProcess14
  - gdcm::UIDs, [772](#)
- JPEGLosslessNonHierarchicalProcess15Retired
  - gdcm::UIDs, [773](#)
- JPEGLosslessProcess14
  - gdcm::TransferSyntax, [759](#)
- JPEGLosslessProcess14\_1
  - gdcm::TransferSyntax, [759](#)
- JPEGSpectralSelectionHierarchicalProcess2022Retired
  - gdcm::UIDs, [773](#)
- JPEGSpectralSelectionHierarchicalProcess2123Retired
  - gdcm::UIDs, [773](#)
- JPEGSpectralSelectionNonHierarchicalProcess68Retired
  - gdcm::UIDs, [772](#)
- JPEGSpectralSelectionNonHierarchicalProcess79Retired
  - gdcm::UIDs, [772](#)
- JPEGSpectralSelectionProcess6\_8
  - gdcm::TransferSyntax, [759](#)

- JPIPIReferenced
  - gdcm::TransferSyntax, [760](#)
  - gdcm::UIDs, [773](#)
- JPIPIReferencedDeflate
  - gdcm::UIDs, [773](#)
- JSON
  - gdcm::JSON, [491](#)
- Join
  - gdcm::Filename, [385](#)
- JunkAfterDocElementError
  - gdcm::Parser, [555](#)
- KAKADUCodec
  - gdcm::KAKADUCodec, [493](#)
- KeyField
  - gdcm::CSAElement, [268](#)
- KeyObjectSelectionDocument
  - gdcm::MediaStorage, [509](#)
- KeyObjectSelectionDocumentStorage
  - gdcm::UIDs, [777](#)
- KeyValuePairArrayType
  - gdcm::CompositeNetworkFunctions, [254](#)
- KeyValuePairType
  - gdcm::CompositeNetworkFunctions, [254](#)
- LD\_ALL
  - gdcm, [126](#)
- LD\_NOSEQ
  - gdcm, [126](#)
- LD\_NOSHADOW
  - gdcm, [126](#)
- LD\_NOSHADOWSEQ
  - gdcm, [126](#)
- LINE
  - gdcm::MeshPrimitive, [518](#)
- LO
  - gdcm::LO, [495](#)
  - gdcm::VR, [853](#)
- LOADBULKDATA
  - gdcm::XMLPrinter, [925](#)
- LOComp
  - gdcm, [124](#)
- LT
  - gdcm::VR, [853](#)
- LTComp
  - gdcm, [124](#)
- LUT
  - gdcm::Bitmap, [216](#)
  - gdcm::ImageCodec, [435](#)
- LUTPtr
  - gdcm::Bitmap, [212](#)
  - gdcm::ImageCodec, [431](#)
- LeadECGWaveformStorage
  - gdcm::MediaStorage, [509](#)
- Level
  - vtkImageMapToWindowLevelColors2, [904](#)
- ListCharSets
  - gdcm::QueryFactory, [616](#)
- LittleEndian
  - gdcm::SwapCode, [730](#)
- Load
  - gdcm::Directory, [330](#)
- LoadDefault
  - gdcm::CSAHeaderDict, [274](#)
  - gdcm::Dict, [315](#)
  - gdcm::PrivateDict, [605](#)
- LoadDefaults
  - gdcm::Defs, [308](#)
  - gdcm::Dicts, [324](#)
- LoadFromDataElement
  - gdcm::CSAHeader, [272](#)
  - gdcm::PDBHeader, [562](#)
- LoadFromFile
  - gdcm::Defs, [308](#)
- LoadIconImage
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMImageReader2, [870](#)
- LoadImageFromFiles
  - gdcm::DirectoryHelper, [332](#)
- LoadOverlays
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMImageReader2, [870](#)
- LoadResourcesFiles
  - gdcm::Global, [406](#)
- LoadSingleFile
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [868](#)
- Locate
  - gdcm::Global, [406](#)
- LodModeType
  - gdcm, [126](#)
- LookupTable
  - gdcm::LookupTable, [498](#)
  - vtkImageMapToColors16, [902](#)
- LookupTableType
  - gdcm::LookupTable, [498](#)
- LossyFlag
  - gdcm::Bitmap, [216](#)
  - gdcm::ImageCodec, [435](#)
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMImageReader2, [870](#)
- m\_ConstMemberFunction
  - gdcm::MemberCommand, [515](#)
- m\_MemberFunction
  - gdcm::MemberCommand, [516](#)
  - gdcm::SimpleMemberCommand, [683](#)
- m\_This
  - gdcm::MemberCommand, [516](#)

- gdcm::SimpleMemberCommand, 683
- m\_char
  - gdcm::ignore\_char, 413
- MAGNIFIED
  - gdcm::Spacing, 695
- MANUAL
  - gdcm::Segment, 650
- mAction
  - gdcm::network::Transition, 763
- MD5
  - gdcm::MD5, 505
- MD5DataImagesType
  - gdcm::Testing, 750
- MD5MetalImagesType
  - vtkGDCMTesting, 885
- mDataSet
  - gdcm::BaseRootQuery, 204
- mElementOffsets
  - gdcm::StreamImageWriter, 706
- mElementOffsets1
  - gdcm::StreamImageWriter, 706
- mEnd
  - gdcm::network::Transition, 763
- mHelpDescription
  - gdcm::BaseRootQuery, 204
- mImage
  - gdcm::BaseRootQuery, 204
- mImplicit
  - gdcm::network::ULConnectionCallback, 822
- MONOCHROME1
  - gdcm::PhotometricInterpretation, 570
- MONOCHROME2
  - gdcm::PhotometricInterpretation, 570
- MPEG2MainProfile
  - gdcm::TransferSyntax, 760
- MPEG2MainProfileMainLevel
  - gdcm::UIDs, 773
- MPTType
  - gdcm::MeshPrimitive, 518
- MPTType\_END
  - gdcm::MeshPrimitive, 518
- mPatient
  - gdcm::BaseRootQuery, 204
- MRImageStorage
  - gdcm::MediaStorage, 508
  - gdcm::UIDs, 775
- MRSpectroscopyStorage
  - gdcm::MediaStorage, 508
  - gdcm::UIDs, 775
- mRootType
  - gdcm::BaseRootQuery, 204
- MS\_END
  - gdcm::MediaStorage, 510
- MSType
  - gdcm::MediaStorage, 508
- mSeries
  - gdcm::BaseRootQuery, 204
- mStudy
  - gdcm::BaseRootQuery, 204
- mWriter
  - gdcm::StreamImageWriter, 707
- mXMax
  - gdcm::StreamImageWriter, 707
- mXMin
  - gdcm::StreamImageWriter, 707
- mYMax
  - gdcm::StreamImageWriter, 707
- mYMin
  - gdcm::StreamImageWriter, 707
- mZMax
  - gdcm::StreamImageWriter, 707
- mZMin
  - gdcm::StreamImageWriter, 707
- Macro
  - gdcm::Macro, 501
- MacroEntry
  - gdcm, 124
- Macros
  - gdcm::Macros, 503
- magenta
  - gdcm::terminal, 138
- MakeDirectory
  - gdcm::System, 735
- MakeNew
  - gdcm::network::Transition, 763
- MakeObject
  - gdcm::AnonymizeEvent, 153
  - gdcm::DataEvent, 294
  - gdcm::DataSetEvent, 303
  - gdcm::Event, 358
  - gdcm::FileNameEvent, 388
  - gdcm::ProgressEvent, 609
- MammographyCADSR
  - gdcm::MediaStorage, 509
- MammographyCADSRStorage
  - gdcm::UIDs, 776
- Mandatory
  - gdcm::Usage, 835
- MapCSAHeaderDictEntry
  - gdcm::CSAHeaderDict, 273
- MapDictEntry
  - gdcm::Dict, 314
- MapIODEntry
  - gdcm::IOD, 459
- MapModuleEntry
  - gdcm::Macro, 501
  - gdcm::Module, 521
- MapScalarsThroughTable2

- vtkLookupTable16, [911](#)
- MapTableEntry
  - gdcm::Table, [737](#)
- MappingType
  - gdcm::Scanner, [645](#)
- MaxLength
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::PersonName, [567](#)
- MaxNumberOfComponents
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::PersonName, [567](#)
- MaxPrintLength
  - gdcm::Printer, [604](#)
- MaximumLengthSub
  - gdcm::network::MaximumLengthSub, [504](#)
- MediaCreationManagementSOPClassUID
  - gdcm::UIDs, [775](#)
- MediaStorage
  - gdcm::MediaStorage, [511](#)
- MediaStorageDataFilesType
  - gdcm::Testing, [750](#)
- MediaStorageDirectoryStorage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [773](#)
- MedicalImageProperties
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMPolyDataReader, [880](#)
  - vtkGDCMPolyDataWriter, [883](#)
- MemberCommand
  - gdcm::MemberCommand, [515](#)
- MeshPrimitive
  - gdcm::MeshPrimitive, [519](#)
- MessageID
  - gdcm::network::CEchoRQ, [234](#)
- MetaInformationTS
  - gdcm::FileMetaInformation, [384](#)
- ModalityPerformedProcedureStepNotificationSOPClass
  - gdcm::UIDs, [774](#)
- ModalityPerformedProcedureStepRetrieveSOPClass
  - gdcm::UIDs, [774](#)
- ModalityPerformedProcedureStepSOPClass
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [774](#)
- ModalityWorklistInformationModelFIND
  - gdcm::UIDs, [777](#)
- Mode
  - gdcm::terminal, [138](#)
- Module
  - gdcm::Module, [522](#)
- ModuleEntry
  - gdcm::ModuleEntry, [524](#)
- ModuleMapType
  - gdcm::Macros, [503](#)
  - gdcm::Modules, [526](#)
- Modules
  - gdcm::Modules, [526](#)
- MovePatientRootQuery
  - gdcm::MovePatientRootQuery, [528](#)
- MoveStudyRootQuery
  - gdcm::MoveStudyRootQuery, [530](#)
- mSPFile
  - gdcm::StreamImageWriter, [707](#)
- MultiframeGrayscaleByteSecondaryCaptureImage↔Storage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- MultiframeGrayscaleWordSecondaryCaptureImage↔Storage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- MultiframeSingleBitSecondaryCaptureImageStorage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- MultiframeTrueColorSecondaryCaptureImageStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [775](#)
- N\_ACTION\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_ACTION\_RSP
  - gdcm::network::DIMSE, [326](#)
- N\_CREATE\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_CREATE\_RSP
  - gdcm::network::DIMSE, [326](#)
- N\_DELETE\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_DELETE\_RSP
  - gdcm::network::DIMSE, [326](#)
- N\_EVENT\_REPORT\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_EVENT\_REPORT\_RSP
  - gdcm::network::DIMSE, [326](#)
- N\_GET\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_GET\_RSP
  - gdcm::network::DIMSE, [326](#)
- N\_SET\_RQ
  - gdcm::network::DIMSE, [326](#)
- N\_SET\_RSP
  - gdcm::network::DIMSE, [326](#)
- NO
  - gdcm::Surface, [719](#)
- NO\_COMPRESSION
  - vtkGDCMImageWriter, [873](#)
- NOMAGIC
  - gdcm::CSAHeader, [270](#)
- Name



- gdcmm::ModuleEntry, 525
- NameField
  - gdcmm::CSAElement, 268
  - gdcmm::PDSElement, 560
- NeedByteSwap
  - gdcmm::Bitmap, 216
  - gdcmm::ImageCodec, 435
- NeedOverlayCleanup
  - gdcmm::ImageCodec, 435
- NegotiatedType
  - gdcmm::TransferSyntax, 759
- NestedMacroEntries
  - gdcmm, 124
- NestedModuleEntries
  - gdcmm::NestedModuleEntries, 533
- New
  - gdcmm::Anonymizer, 157
  - gdcmm::FileChangeTransferSyntax, 374
  - gdcmm::FileStreamer, 394
  - gdcmm::MemberCommand, 515
  - gdcmm::Scanner, 647
  - gdcmm::SequenceOfFragments, 663
  - gdcmm::SequenceOfItems, 668
  - gdcmm::ServiceClassUser, 677
  - gdcmm::SimpleMemberCommand, 683
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
  - vtkGDCMImageWriter, 873
  - vtkGDCMMedicalImageProperties, 877
  - vtkGDCMPolyDataReader, 879
  - vtkGDCMPolyDataWriter, 882
  - vtkGDCMTesting, 885
  - vtkGDCMThreadedImageReader, 888
  - vtkGDCMThreadedImageReader2, 890
  - vtkImageColorViewer, 896
  - vtkImageMapToColors16, 901
  - vtkImageMapToWindowLevelColors2, 904
  - vtkImagePlanarComponentsToComponents, 906
  - vtkImageRGBToYBR, 907
  - vtkImageYBRToRGB, 909
  - vtkLookupTable16, 911
  - vtkRTStructSetProperties, 914
- NoElementsError
  - gdcmm::Parser, 555
- NoError
  - gdcmm::Parser, 555
- NoMemoryError
  - gdcmm::Parser, 555
- NoObject
  - gdcmm::MediaStorage, 510
- NoOfItemsField
  - gdcmm::CSAElement, 268
- Normalize
  - gdcmm::DirectionCosines, 327
- NuclearMedicineImageStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 776
- NuclearMedicineImageStorageRetired
  - gdcmm::MediaStorage, 508
  - gdcmm::UIDs, 775
- NumberOfDimensions
  - gdcmm::Bitmap, 216
  - gdcmm::ImageCodec, 435
- NumberOfIconImages
  - vtkGDCMImageReader, 865
  - vtkGDCMImageReader2, 870
- NumberOfOverlays
  - vtkGDCMImageReader, 865
  - vtkGDCMImageReader2, 870
- NumberOfSurfaces
  - gdcmm::SurfaceWriter, 729
- OB
  - gdcmm::VR, 853
- OB\_OW
  - gdcmm::VR, 853
- OBLIQUE
  - gdcmm::Orientation, 545
- OD
  - gdcmm::VR, 853
- OF
  - gdcmm::VR, 853
- OPENSSL
  - gdcmm::CryptoFactory, 261
- OPENSSL7
  - gdcmm::CryptoFactory, 261
- OW
  - gdcmm::VR, 853
- Object
  - gdcmm::Object, 536
- ObjectEnd
  - gdcmm::MediaStorage, 510
- ObjectType
  - gdcmm::MediaStorage, 510
- Ofstream
  - gdcmm::Writer, 921
- OnlyUUID
  - gdcmm::XMLPrinter, 925
- op
  - gdcmm::SerieHelper::Rule, 642
- OpenSSLCryptoFactory
  - gdcmm::OpenSSLCryptoFactory, 538
- OpenSSLCryptographicMessageSyntax
  - gdcmm::OpenSSLCryptographicMessageSyntax, 539
- OpenSSL7CryptoFactory
  - gdcmm::OpenSSL7CryptoFactory, 541
- OpenSSL7CryptographicMessageSyntax



- gdcmm::OpenSSLP7CryptographicMessageSyntax, 543
- operator const char \*
  - gdcmm::ConstCharWrapper, 258
  - gdcmm::Filename, 386
  - gdcmm::String, 710
- operator const double \*
  - gdcmm::DirectionCosines, 328
- operator const std::vector< char > &
  - gdcmm::ByteValue, 228
- operator MStype
  - gdcmm::MediaStorage, 512
- operator ObjectType \*
  - gdcmm::SmartPointer, 687
- operator PIType
  - gdcmm::PhotometricInterpretation, 571
- operator ScalarType
  - gdcmm::PixelFormat, 575
- operator SwapCode::SwapCodeType
  - gdcmm::SwapCode, 731
- operator TSType
  - gdcmm::TransferSyntax, 761
  - gdcmm::UIDs, 786
- operator TypeType
  - gdcmm::Type, 765
- operator uint32\_t
  - gdcmm::VL, 846
- operator UsageType
  - gdcmm::Usage, 835
- operator VMType
  - gdcmm::VM, 850
- operator VRType
  - gdcmm::VR, 854
- operator!=
  - gdcmm, 126
  - gdcmm::Attribute, 173
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 179
  - gdcmm::CodeString, 248
  - gdcmm::PixelFormat, 575
  - gdcmm::Tag, 745
- operator<
  - gdcmm::Attribute, 173
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 179
  - gdcmm::CSAElement, 267
  - gdcmm::CSAHeaderDictEntry, 275
  - gdcmm::DataElement, 288
  - gdcmm::PrivateTag, 607
  - gdcmm::Tag, 745
- operator<<
  - gdcmm, 126–130
  - gdcmm::BasicOffsetTable, 209
  - gdcmm::CSAElement, 268
- gdcmm::CSAHeader, 272
- gdcmm::CSAHeaderDict, 274
- gdcmm::CSAHeaderDictEntry, 276
- gdcmm::CodeString, 248
- gdcmm::CommandDataSet, 252
- gdcmm::DataElement, 291
- gdcmm::DataSet, 301
- gdcmm::Dict, 315
- gdcmm::DictEntry, 320
- gdcmm::Dicts, 324
- gdcmm::Directory, 330
- gdcmm::File, 369
- gdcmm::FileMetaInformation, 383
- gdcmm::FileSet, 392
- gdcmm::Fragment, 404
- gdcmm::Global, 406
- gdcmm::GroupDict, 408
- gdcmm::IOD, 460
- gdcmm::IODEntry, 462
- gdcmm::IODs, 463
- gdcmm::Item, 470
- gdcmm::Macro, 502
- gdcmm::Macros, 503
- gdcmm::MediaStorage, 512
- gdcmm::Module, 522
- gdcmm::ModuleEntry, 525
- gdcmm::Modules, 526
- gdcmm::NestedModuleEntries, 533
- gdcmm::Object, 536
- gdcmm::Orientation, 546
- gdcmm::PDBelement, 560
- gdcmm::PDBHeader, 562
- gdcmm::PhotometricInterpretation, 571
- gdcmm::PixelFormat, 576
- gdcmm::Preamble, 592
- gdcmm::PrivateDict, 605
- gdcmm::PrivateTag, 607
- gdcmm::Scanner, 648
- gdcmm::Sorter, 694
- gdcmm::SwapCode, 731
- gdcmm::Table, 737
- gdcmm::Tag, 748
- gdcmm::TransferSyntax, 761
- gdcmm::Type, 765
- gdcmm::UI, 766
- gdcmm::Usage, 835
- gdcmm::VL, 846
- gdcmm::VM, 850
- gdcmm::VR, 855
- gdcmm::Version, 844
- operator<=
  - gdcmm::Tag, 745
- operator>>
  - gdcmm, 130

- gdcmm::Tag, 748
- operator\*
  - gdcmm::SmartPointer, 688
- operator()
  - gdcmm::DataSet, 300
  - gdcmm::Scanner::ltstr, 500
- operator++
  - gdcmm::VL, 846
- operator+=
  - gdcmm::VL, 846
- operator->
  - gdcmm::SmartPointer, 688
- operator=
  - gdcmm::BoxRegion, 221
  - gdcmm::ByteValue, 228
  - gdcmm::CSAElement, 267
  - gdcmm::DataElement, 288
  - gdcmm::DataSet, 300
  - gdcmm::Element< TVR, VM::VM1\_n >, 341
  - gdcmm::Object, 536
  - gdcmm::ParseException, 553
  - gdcmm::Preamble, 592
  - gdcmm::SequenceOfItems, 668
  - gdcmm::SmartPointer, 688
  - gdcmm::Tag, 745
  - gdcmm::network::UserInformation, 838
- operator==
  - gdcmm, 130
  - gdcmm::Attribute, 173
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 179
  - gdcmm::ByteValue, 228
  - gdcmm::CSAElement, 267
  - gdcmm::CodeString, 248
  - gdcmm::DataElement, 288
  - gdcmm::PDBelement, 560
  - gdcmm::PixelFormat, 575, 576
  - gdcmm::PresentationContext, 594
  - gdcmm::SequenceOfFragments, 663
  - gdcmm::SequenceOfItems, 668
  - gdcmm::Tag, 745
  - gdcmm::Value, 842
  - gdcmm::network::AbstractSyntax, 151
  - gdcmm::network::PresentationContextRQ, 599
  - gdcmm::network::TransferSyntaxSub, 762
- operator[]
  - gdcmm::Attribute, 173, 174
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, 185
  - gdcmm::DataSet, 300
  - gdcmm::Element, 337
  - gdcmm::Element< TVR, VM::VM1\_n >, 341
  - gdcmm::Tag, 746
- OphthalmicPhotography16BitImageStorage
  - gdcmm::UIDs, 776
- OphthalmicPhotography8BitImageStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 776
- OphthalmicTomographyImageStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 776
- OrderFileList
  - gdcmm::SerieHelper, 672
- Orientation
  - gdcmm::Orientation, 545
- OrientationType
  - gdcmm::Orientation, 545
- Output
  - gdcmm::BitmapToBitmapFilter, 218
- OutputFormat
  - vtkImageMapToColors16, 902
- OutputTypes
  - gdcmm::DictConverter, 317
- Overlay
  - gdcmm::Overlay, 549
- OverlayImageActor
  - vtkImageColorViewer, 899
- OverlayType
  - gdcmm::Overlay, 549
- Overlays
  - gdcmm::Pixmap, 579
- PALETTE\_COLOR
  - gdcmm::PhotometricInterpretation, 570
- PDBelement
  - gdcmm::PDBelement, 560
- PDBHeader
  - gdcmm::PDBHeader, 562
- PDF
  - gdcmm::MediaStorage, 510
- PDFCodec
  - gdcmm::PDFCodec, 564
- PDataTFPDU
  - gdcmm::network::PDataTFPDU, 558
- PETImageStorage
  - gdcmm::MediaStorage, 509
- PF
  - gdcmm::Bitmap, 216
  - gdcmm::ImageCodec, 435
- PGXCodec
  - gdcmm::PGXCodec, 568
- PHILIPS
  - gdcmm::Dicts, 323
- PI
  - gdcmm::Bitmap, 216
  - gdcmm::ImageCodec, 435
- PI\_END
  - gdcmm::PhotometricInterpretation, 570

- PType
  - gdcm::PhotometricInterpretation, [570](#)
- PN
  - gdcm::VR, [853](#)
- PNComp
  - gdcm, [124](#)
- PNMCodec
  - gdcm::PNMCodec, [590](#)
- POINTS
  - gdcm::Surface, [719](#)
- PVRGCodec
  - gdcm::PVRGCodec, [611](#)
- Pack
  - gdcm::Unpacker12Bits, [834](#)
- Padding
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::PersonName, [567](#)
- Parent
  - gdcm::Element< TVR, VM::VM1\_2 >, [339](#)
  - gdcm::Element< TVR, VM::VM2\_2n >, [343](#)
  - gdcm::Element< TVR, VM::VM2\_n >, [344](#)
  - gdcm::Element< TVR, VM::VM3\_3n >, [346](#)
  - gdcm::Element< TVR, VM::VM3\_n >, [347](#)
- Parse
  - gdcm::Parser, [555](#)
- ParseBuffer
  - gdcm::Parser, [556](#)
- ParseCertificateFile
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [540](#)
  - gdcm::OpenSSL7CryptographicMessageSyntax, [543](#)
- ParseDateTime
  - gdcm::System, [735](#), [736](#)
- ParseDump
  - gdcm::ASN1, [168](#)
- ParseDumpFile
  - gdcm::ASN1, [168](#)
- ParseException
  - gdcm::ParseException, [553](#)
- ParseKeyFile
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [540](#)
  - gdcm::OpenSSL7CryptographicMessageSyntax, [544](#)
- Parser
  - gdcm::Parser, [555](#)
- PassAlphaToOutput
  - vtkImageMapToColors16, [902](#)
- Patient
  - gdcm::Patient, [556](#)
- PatientRootQueryRetrieveInformationModelFIND
  - gdcm::UIDs, [777](#)
- PatientRootQueryRetrieveInformationModelGET
  - gdcm::UIDs, [777](#)
- PatientRootQueryRetrieveInformationModelMOVE
  - gdcm::UIDs, [777](#)
- PatientStudyOnlyQueryRetrieveInformationModelFIND↔
  - Retired
  - gdcm::UIDs, [777](#)
- PatientStudyOnlyQueryRetrieveInformationModelGET↔
  - Retired
  - gdcm::UIDs, [777](#)
- PatientStudyOnlyQueryRetrieveInformationModelMOV↔
  - ERetired
  - gdcm::UIDs, [777](#)
- PerformAction
  - gdcm::network::ULAction, [788](#)
  - gdcm::network::ULActionAA1, [789](#)
  - gdcm::network::ULActionAA2, [790](#)
  - gdcm::network::ULActionAA3, [791](#)
  - gdcm::network::ULActionAA4, [792](#)
  - gdcm::network::ULActionAA5, [793](#)
  - gdcm::network::ULActionAA6, [794](#)
  - gdcm::network::ULActionAA7, [795](#)
  - gdcm::network::ULActionAA8, [796](#)
  - gdcm::network::ULActionAE1, [797](#)
  - gdcm::network::ULActionAE2, [798](#)
  - gdcm::network::ULActionAE3, [799](#)
  - gdcm::network::ULActionAE4, [800](#)
  - gdcm::network::ULActionAE5, [801](#)
  - gdcm::network::ULActionAE6, [802](#)
  - gdcm::network::ULActionAE7, [803](#)
  - gdcm::network::ULActionAE8, [804](#)
  - gdcm::network::ULActionAR1, [805](#)
  - gdcm::network::ULActionAR10, [806](#)
  - gdcm::network::ULActionAR2, [807](#)
  - gdcm::network::ULActionAR3, [808](#)
  - gdcm::network::ULActionAR4, [809](#)
  - gdcm::network::ULActionAR5, [810](#)
  - gdcm::network::ULActionAR6, [811](#)
  - gdcm::network::ULActionAR7, [812](#)
  - gdcm::network::ULActionAR8, [813](#)
  - gdcm::network::ULActionAR9, [814](#)
  - gdcm::network::ULActionDT1, [815](#)
  - gdcm::network::ULActionDT2, [816](#)
- Philips3D
  - gdcm::MediaStorage, [509](#)
- PhilipsPrivateMRSyntheticImageStorage
  - gdcm::MediaStorage, [510](#)
- PhotometricInterpretation
  - gdcm::PhotometricInterpretation, [571](#)
- PixelData
  - gdcm::Bitmap, [216](#)
  - gdcm::PixmapReader, [583](#)
  - gdcm::PixmapWriter, [588](#)

- PixelFormat
  - gdcm::PixelFormat, [574](#)
- Pixmap
  - gdcm::Pixmap, [578](#)
- PixmapReader
  - gdcm::Bitmap, [216](#)
  - gdcm::PixmapReader, [582](#)
- PixmapToPixmapFilter
  - gdcm::PixmapToPixmapFilter, [584](#)
- PixmapWriter
  - gdcm::PixmapWriter, [587](#)
- PlanarConfiguration
  - gdcm::Bitmap, [216](#)
  - gdcm::ImageCodec, [435](#)
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMImageReader2, [871](#)
- pointer
  - gdcm::CodeString, [247](#)
  - gdcm::LO, [495](#)
  - gdcm::String, [709](#)
- PositronEmissionTomographyImageStorage
  - gdcm::UIDs, [777](#)
- Preamble
  - gdcm::Preamble, [592](#)
- PrepareWrite
  - gdcm::PixmapWriter, [587](#)
  - gdcm::SegmentWriter, [658](#)
  - gdcm::SurfaceWriter, [729](#)
- PrepareWritePointMacro
  - gdcm::SurfaceWriter, [729](#)
- Prepend
  - gdcm::Global, [406](#)
- PresentationContext
  - gdcm::PresentationContext, [593](#)
- PresentationContextAC
  - gdcm::network::PresentationContextAC, [594](#)
- PresentationContextArrayType
  - gdcm::PresentationContextGenerator, [596](#)
  - gdcm::network::AAssociateRQPDU, [147](#)
- PresentationContextGenerator
  - gdcm::PresentationContextGenerator, [596](#)
- PresentationContextRQ
  - gdcm::network::PresentationContextRQ, [598](#)
- PresentationDataValue
  - gdcm::network::PresentationDataValue, [600](#)
- PresentationLUTSOPClass
  - gdcm::UIDs, [775](#)
- PrettyPrintOff
  - gdcm::JSON, [491](#)
- PrettyPrintOn
  - gdcm::JSON, [491](#)
- PrimitiveData
  - gdcm::MeshPrimitive, [519](#)
- PrimitiveType
  - gdcm::MeshPrimitive, [519](#)
- PrimitivesData
  - gdcm::MeshPrimitive, [518](#)
- Print
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::Attribute, [174](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [179](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [185](#)
  - gdcm::BaseRootQuery, [203](#)
  - gdcm::Bitmap, [214](#)
  - gdcm::BoxRegion, [221](#)
  - gdcm::ByteValue, [228](#)
  - gdcm::CSAHeader, [272](#)
  - gdcm::Curve, [282](#)
  - gdcm::DataSet, [300](#)
  - gdcm::DictPrinter, [322](#)
  - gdcm::DirectionCosines, [328](#)
  - gdcm::Directory, [330](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
  - gdcm::Element< VR::AS, VM::VM5 >, [348](#)
  - gdcm::Event, [358](#)
  - gdcm::Image, [416](#)
  - gdcm::LookupTable, [499](#)
  - gdcm::Object, [536](#)
  - gdcm::Orientation, [546](#)
  - gdcm::Overlay, [551](#)
  - gdcm::PDBHeader, [562](#)
  - gdcm::PersonName, [566](#)
  - gdcm::PixelFormat, [576](#)
  - gdcm::Pixmap, [579](#)
  - gdcm::Preamble, [592](#)
  - gdcm::PresentationContext, [594](#)
  - gdcm::Printer, [603](#)
  - gdcm::Region, [634](#)
  - gdcm::Scanner, [647](#)
  - gdcm::SegmentedPaletteColorLookupTable, [653](#)
  - gdcm::SequenceOfFragments, [663](#)
  - gdcm::SequenceOfItems, [668](#)
  - gdcm::Sorter, [693](#)
  - gdcm::TagPath, [749](#)
  - gdcm::Testing, [753](#)
  - gdcm::Version, [844](#)
  - gdcm::XMLPrinter, [925](#)
  - gdcm::network::AAAbortPDU, [140](#)
  - gdcm::network::AAssociateACPDU, [143](#)
  - gdcm::network::AAssociateRJPDU, [145](#)
  - gdcm::network::AAssociateRQPDU, [148](#)
  - gdcm::network::AReleaseRPPDU, [164](#)
  - gdcm::network::AReleaseRPPDU, [166](#)
  - gdcm::network::AbstractSyntax, [151](#)
  - gdcm::network::ApplicationContext, [161](#)

- gdcm::network::AsynchronousOperationsWindow↔  
Sub, 168
- gdcm::network::BasePDU, 200
- gdcm::network::ImplementationClassUIDSub, 454
- gdcm::network::ImplementationVersionNameSub,  
455
- gdcm::network::MaximumLengthSub, 504
- gdcm::network::PDataTFPDU, 558
- gdcm::network::PresentationContextAC, 595
- gdcm::network::PresentationContextRQ, 599
- gdcm::network::PresentationDataValue, 600
- gdcm::network::RoleSelectionSub, 641
- gdcm::network::SOPClassExtendedNegociationSub,  
689
- gdcm::network::ServiceClassApplicationInformation,  
673
- gdcm::network::TransferSyntaxSub, 762
- gdcm::network::UserInformation, 838
- PrintASCII
  - gdcm::ByteValue, 228
- PrintASCIIXML
  - gdcm::ByteValue, 228
- PrintAsContinuousString
  - gdcm::Tag, 746
- PrintAsContinuousUpperCaseString
  - gdcm::Tag, 746
- PrintAsPipeSeparatedString
  - gdcm::Tag, 746
- PrintDataElement
  - gdcm::Printer, 603
  - gdcm::XMLPrinter, 925
- PrintDataElement2
  - gdcm::DictPrinter, 322
- PrintDataSet
  - gdcm::Printer, 603
  - gdcm::XMLPrinter, 925
- PrintDataSet2
  - gdcm::DictPrinter, 322
- PrintGroupLength
  - gdcm::ByteValue, 228
- PrintHex
  - gdcm::ByteValue, 228
- PrintHexXML
  - gdcm::ByteValue, 228
- PrintJobSOPClass
  - gdcm::UIDs, 774
- PrintPNXML
  - gdcm::ByteValue, 228
- PrintQueueManagementSOPClassRetired
  - gdcm::UIDs, 775
- PrintQueueSOPInstanceRetired
  - gdcm::UIDs, 775
- PrintSQ
  - gdcm::Printer, 603
- gdcm::XMLPrinter, 925
- PrintSelf
  - vtkGDCMImageReader, 862
  - vtkGDCMImageReader2, 868
  - vtkGDCMImageWriter, 873
  - vtkGDCMMedicalImageProperties, 877
  - vtkGDCMPolyDataReader, 879
  - vtkGDCMPolyDataWriter, 882
  - vtkGDCMTesting, 885
  - vtkGDCMThreadedImageReader, 888
  - vtkGDCMThreadedImageReader2, 890
  - vtkImageColorViewer, 896
  - vtkImageMapToColors16, 901
  - vtkImageMapToWindowLevelColors2, 904
  - vtkImagePlanarComponentsToComponents, 906
  - vtkImageRGBToYBR, 907
  - vtkImageYBRToRGB, 909
  - vtkLookupTable16, 911
  - vtkRTStructSetProperties, 915
- PrintStyle
  - gdcm::Printer, 604
  - gdcm::XMLPrinter, 925
- PrintStyles
  - gdcm::Printer, 603
  - gdcm::XMLPrinter, 925
- PrintTable
  - gdcm::network::ULTransitionTable, 828
- PrintXML
  - gdcm::PrivateDict, 605
- Printer
  - gdcm::Printer, 603
- PrinterConfigurationRetrievalSOPClass
  - gdcm::UIDs, 774
- PrinterConfigurationRetrievalSOPInstance
  - gdcm::UIDs, 774
- PrinterSOPClass
  - gdcm::UIDs, 774
- PrinterSOPInstance
  - gdcm::UIDs, 774
- PrivateDict
  - gdcm::PrivateDict, 605
- PrivateTag
  - gdcm::PrivateTag, 607
- ProceduralEventLoggingSOPClass
  - gdcm::UIDs, 774
- ProceduralEventLoggingSOPInstance
  - gdcm::UIDs, 774
- ProcedureLogStorage
  - gdcm::UIDs, 776
- Process
  - gdcm::Parser, 556
- ProcessDataSet
  - gdcm::FileExplicitFilter, 378
- ProcessPublicTag

- gdcmm::Scanner, 647
- ProcessRequest
  - vtkGDCMImageReader2, 868
- ProduceCharacterSetDataElement
  - gdcmm::QueryFactory, 616
- ProduceQuery
  - gdcmm::QueryFactory, 616
- ProductCharacteristicsQuerySOPClass
  - gdcmm::UIDs, 778
- ProgressEvent
  - gdcmm::ProgressEvent, 609
- PropertyCategory
  - gdcmm::Segment, 651
- PropertyType
  - gdcmm::Segment, 651
- PseudoColorSoftcopyPresentationStateStorageSOP↔
  - Class
    - gdcmm::UIDs, 776
- PullPrintRequestSOPClassRetired
  - gdcmm::UIDs, 775
- PullStoredPrintManagementMetaSOPClassRetired
  - gdcmm::UIDs, 775
- Push
  - gdcmm::TagPath, 749
- PushBackFile
  - vtkGDCMMedicalImageProperties, 877
- PythonFilter
  - gdcmm::PythonFilter, 613
- Quality
  - gdcmm::JPEGCodec, 486
- QueryFactory
  - gdcmm::BaseRootQuery, 204
  - gdcmm::FindPatientRootQuery, 399
  - gdcmm::FindStudyRootQuery, 401
  - gdcmm::MovePatientRootQuery, 529
  - gdcmm::MoveStudyRootQuery, 531
- RAWCodec
  - gdcmm::RAWCodec, 626
- README.txt, 1192
- RED
  - gdcmm::LookupTable, 498
- RFC2557MIMEencapsulation
  - gdcmm::UIDs, 773
- RGB
  - gdcmm::PhotometricInterpretation, 570
- RGB2YBR
  - gdcmm::ImageChangePhotometricInterpretation, 421
- RGBPixelsToRGBPlanes
  - gdcmm::ImageChangePlanarConfiguration, 424
- RGBPlanesToRGBPixels
  - gdcmm::ImageChangePlanarConfiguration, 424
- RGBToRecommendedDisplayCIELab
  - gdcmm::SurfaceHelper, 724
- RGBToRecommendedDisplayGrayscale
  - gdcmm::SurfaceHelper, 724
- RLE\_COMPRESSION
  - vtkGDCMImageWriter, 873
- RLECodec
  - gdcmm::RLECodec, 638
- RLELossless
  - gdcmm::TransferSyntax, 760
  - gdcmm::UIDs, 773
- ROI
  - gdcmm::Overlay, 549
- RTBeamsDeliveryInstructionStorageSupplement74↔
  - FrozenDraft
    - gdcmm::UIDs, 777
- RTBeamsTreatmentRecordStorage
  - gdcmm::UIDs, 777
- RTBrachyTreatmentRecordStorage
  - gdcmm::UIDs, 777
- RTConventionalMachineVerificationSupplement74↔
  - FrozenDraft
    - gdcmm::UIDs, 777
- RTDoseStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 777
- RTImageStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 777
- RTIonBeamsTreatmentRecordStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 777
- RTIonMachineVerificationSupplement74FrozenDraft
  - gdcmm::UIDs, 777
- RTIonPlanStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 777
- RTPlanStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 777
- RTStructSetProperties
  - vtkGDCMPolyDataReader, 880
  - vtkGDCMPolyDataWriter, 883
- RTStructureSetStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 777
- RTTreatmentSummaryRecordStorage
  - gdcmm::MediaStorage, 510
  - gdcmm::UIDs, 777
- RawDataStorage
  - gdcmm::MediaStorage, 509
  - gdcmm::UIDs, 776
- Read
  - gdcmm::BasicOffsetTable, 208
  - gdcmm::ByteValue, 228, 229
  - gdcmm::CP246ExplicitDataElement, 259



- gdcm::CSAHeader, [272](#)
- gdcm::CommandDataSet, [252](#)
- gdcm::DataElement, [288](#)
- gdcm::DataSet, [300](#)
- gdcm::Element, [337](#)
- gdcm::Element< TVR, VM::VM1\_n >, [341](#)
- gdcm::EncodingImplementation< VR::VRASCII >, [353](#)
- gdcm::EncodingImplementation< VR::VRBINARY >, [354](#)
- gdcm::ExplicitDataElement, [363](#)
- gdcm::ExplicitImplicitDataElement, [365](#)
- gdcm::File, [368](#)
- gdcm::FileMetaInformation, [382](#)
- gdcm::Fragment, [403](#)
- gdcm::ImageReader, [445](#)
- gdcm::ImageRegionReader, [448](#)
- gdcm::ImplicitDataElement, [457](#)
- gdcm::Item, [470](#)
- gdcm::PGXCodec, [569](#)
- gdcm::PNMCodec, [590](#)
- gdcm::PixmapReader, [582](#)
- gdcm::Preamble, [592](#)
- gdcm::Reader, [630](#)
- gdcm::SegmentReader, [656](#)
- gdcm::SequenceOfFragments, [663](#)
- gdcm::SequenceOfItems, [669](#)
- gdcm::StreamImageReader, [701](#)
- gdcm::SurfaceReader, [727](#)
- gdcm::TableReader, [740](#)
- gdcm::Tag, [746](#)
- gdcm::UNExplicitDataElement, [831](#)
- gdcm::UNExplicitImplicitDataElement, [833](#)
- gdcm::VL, [846](#)
- gdcm::VR, [855](#)
- gdcm::VR16ExplicitDataElement, [856](#)
- gdcm::VRVLSize< 0 >, [858](#)
- gdcm::VRVLSize< 1 >, [858](#)
- gdcm::ValueIO, [843](#)
- gdcm::network::AAAbortPDU, [140](#)
- gdcm::network::AAssociateACPDU, [143](#)
- gdcm::network::AAssociateRJPDU, [145](#)
- gdcm::network::AAssociateRQPDU, [148](#)
- gdcm::network::AReleaseRPPDU, [164](#)
- gdcm::network::AReleaseRQPDU, [166](#)
- gdcm::network::AbstractSyntax, [151](#)
- gdcm::network::ApplicationContext, [161](#)
- gdcm::network::AsynchronousOperationsWindow↵Sub, [169](#)
- gdcm::network::BasePDU, [200](#)
- gdcm::network::ImplementationClassUIDSub, [454](#)
- gdcm::network::ImplementationVersionNameSub, [455](#)
- gdcm::network::MaximumLengthSub, [504](#)
- gdcm::network::PDataTFPDU, [558](#)
- gdcm::network::PresentationContextAC, [595](#)
- gdcm::network::PresentationContextRQ, [599](#)
- gdcm::network::PresentationDataValue, [600](#)
- gdcm::network::RoleSelectionSub, [641](#)
- gdcm::network::SOPClassExtendedNegociationSub, [689](#)
- gdcm::network::ServiceClassApplicationInformation, [673](#)
- gdcm::network::TransferSyntaxSub, [762](#)
- gdcm::network::UserInformation, [838](#)
- Read16
  - gdcm::VL, [846](#)
- ReadACRNEMAIImage
  - gdcm::ImageReader, [446](#)
  - gdcm::PixmapReader, [582](#)
- ReadBacktrack
  - gdcm::Fragment, [403](#)
- ReadCompat
  - gdcm::FileMetaInformation, [382](#)
- ReadCompatInternal
  - gdcm::FileMetaInformation, [383](#)
- ReadComputeLength
  - gdcm::EncodingImplementation< VR::VRASCII >, [353](#)
  - gdcm::EncodingImplementation< VR::VRBINARY >, [354](#)
- ReadDataSet
  - gdcm::Reader, [631](#)
- ReadFiles
  - vtkGDCMThreadedImageReader, [888](#)
- ReadFromCommaSeparatedString
  - gdcm::PrivateTag, [607](#)
  - gdcm::Tag, [746](#)
- ReadFromContinuousString
  - gdcm::Tag, [746](#)
- ReadFromPipeSeparatedString
  - gdcm::Tag, [746](#)
- ReadImage
  - gdcm::ImageReader, [446](#)
  - gdcm::PixmapReader, [582](#)
- ReadImageInformation
  - gdcm::StreamImageReader, [701](#)
- ReadImageInternal
  - gdcm::PixmapReader, [582](#)
- ReadInformation
  - gdcm::ImageRegionReader, [448](#)
- ReadInto
  - gdcm::network::PDataTFPDU, [558](#)
  - gdcm::network::PresentationDataValue, [600](#)
- ReadIntoBuffer
  - gdcm::ImageRegionReader, [448](#)
- ReadMetaInformation
  - gdcm::Reader, [631](#)

- ReadNested
  - gdcm::DataSet, [300](#)
- ReadNoSwap
  - gdcm::EncodingImplementation< VR::VRASCII >, [354](#)
  - gdcm::EncodingImplementation< VR::VRBINARY >, [355](#)
- ReadOrSkip
  - gdcm::DataElement, [289](#)
- ReadPointMacro
  - gdcm::SurfaceReader, [727](#)
- ReadPreValue
  - gdcm::CP246ExplicitDataElement, [259](#)
  - gdcm::DataElement, [289](#)
  - gdcm::ExplicitDataElement, [363](#)
  - gdcm::ExplicitImplicitDataElement, [365](#)
  - gdcm::Fragment, [403](#)
  - gdcm::ImplicitDataElement, [457](#)
  - gdcm::SequenceOfFragments, [663](#)
  - gdcm::UNExplicitDataElement, [831](#)
  - gdcm::UNExplicitImplicitDataElement, [833](#)
  - gdcm::VR16ExplicitDataElement, [856](#)
- ReadPreamble
  - gdcm::Reader, [631](#)
- ReadSegment
  - gdcm::SegmentReader, [656](#)
- ReadSegments
  - gdcm::SegmentReader, [656](#)
- ReadSelectedPrivateTags
  - gdcm::DataSet, [300](#)
  - gdcm::Reader, [631](#)
- ReadSelectedPrivateTagsWithLength
  - gdcm::DataSet, [300](#)
- ReadSelectedTags
  - gdcm::DataSet, [300](#)
  - gdcm::Reader, [631](#)
- ReadSelectedTagsWithLength
  - gdcm::DataSet, [300](#)
- ReadSurface
  - gdcm::SurfaceReader, [727](#)
- ReadSurfaces
  - gdcm::SurfaceReader, [727](#)
- ReadUpToTag
  - gdcm::DataSet, [300](#)
  - gdcm::Reader, [631](#)
- ReadUpToTagWithLength
  - gdcm::DataSet, [300](#)
- ReadVM
  - gdcm::DictConverter, [317](#)
- ReadVR
  - gdcm::DictConverter, [317](#)
- ReadValue
  - gdcm::CP246ExplicitDataElement, [260](#)
  - gdcm::DataElement, [289](#)
  - gdcm::ExplicitDataElement, [363](#)
  - gdcm::ExplicitImplicitDataElement, [365](#)
  - gdcm::Fragment, [403](#)
  - gdcm::ImplicitDataElement, [457](#)
  - gdcm::SequenceOfFragments, [663](#)
  - gdcm::UNExplicitDataElement, [831](#)
  - gdcm::UNExplicitImplicitDataElement, [833](#)
  - gdcm::VR16ExplicitDataElement, [857](#)
- ReadValueWithLength
  - gdcm::DataElement, [289](#)
  - gdcm::ImplicitDataElement, [457](#)
- ReadWithLength
  - gdcm::CP246ExplicitDataElement, [260](#)
  - gdcm::DataElement, [289](#)
  - gdcm::DataSet, [300](#)
  - gdcm::ExplicitDataElement, [363](#)
  - gdcm::ExplicitImplicitDataElement, [365](#)
  - gdcm::ImplicitDataElement, [457](#)
  - gdcm::UNExplicitDataElement, [831](#)
  - gdcm::VR16ExplicitDataElement, [857](#)
- Reader
  - gdcm::Reader, [630](#)
- Readuint16
  - gdcm::DictConverter, [317](#)
- RealWorldValueMappingStorage
  - gdcm::UIDs, [776](#)
- RecommendedDisplayCIELabToRGB
  - gdcm::SurfaceHelper, [723](#)
- RecurseDataSet
  - gdcm::Anonymizer, [157](#)
- red
  - gdcm::terminal, [137](#)
- reference
  - gdcm::CodeString, [247](#)
  - gdcm::LO, [495](#)
  - gdcm::String, [709](#)
- ReferenceFrameOfReferenceUID
  - vtkRTStructSetProperties, [915](#)
- ReferenceSeriesInstanceUID
  - vtkRTStructSetProperties, [915](#)
- ReferencedColorPrintManagementMetaSOPClassRetired
  - gdcm::UIDs, [774](#)
- ReferencedGrayscalePrintManagementMetaSOPClass↔Retired
  - gdcm::UIDs, [774](#)
- ReferencedImageBoxSOPClassRetired
  - gdcm::UIDs, [774](#)
- Region
  - gdcm::Region, [633](#)
- Register
  - gdcm::Object, [536](#)
- Remove
  - gdcm::Anonymizer, [157](#)
  - gdcm::DataSet, [300](#)



- gdcm::FileAnonymizer, [371](#)
- gdcm::Preamble, [592](#)
- RemoveAllObservers
  - gdcm::Subject, [715](#)
- RemoveDictEntry
  - gdcm::PrivateDict, [605](#)
- RemoveFile
  - gdcm::System, [736](#)
- RemoveGroupLength
  - gdcm::Anonymizer, [157](#)
- RemoveItemByIndex
  - gdcm::SequenceOfItems, [669](#)
- RemoveObserver
  - gdcm::Subject, [716](#)
- RemoveOverlay
  - gdcm::Pixmap, [579](#)
- RemovePrivateTags
  - gdcm::Anonymizer, [157](#)
- RemoveRetired
  - gdcm::Anonymizer, [158](#)
- Render
  - vtkImageColorViewer, [896](#)
- RenderWindow
  - vtkImageColorViewer, [899](#)
- Renderer
  - vtkImageColorViewer, [899](#)
- Replace
  - gdcm::Anonymizer, [158](#)
  - gdcm::CommandDataSet, [252](#)
  - gdcm::DataSet, [300](#)
  - gdcm::FileAnonymizer, [371](#)
  - gdcm::FileMetaInformation, [383](#)
- ReplaceEmpty
  - gdcm::DataSet, [301](#)
- RequestData
  - vtkGDCMImageReader2, [868](#)
  - vtkGDCMPolyDataReader, [879](#)
  - vtkImageMapToColors16, [901](#)
  - vtkImageMapToWindowLevelColors2, [904](#)
  - vtkImagePlanarComponentsToComponents, [906](#)
- RequestData\_HemodynamicWaveformStorage
  - vtkGDCMPolyDataReader, [879](#)
- RequestData\_RTStructureSetStorage
  - vtkGDCMPolyDataReader, [879](#)
- RequestDataCompat
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [868](#)
  - vtkGDCMThreadedImageReader, [888](#)
- RequestInformation
  - vtkGDCMImageReader2, [868](#)
  - vtkGDCMPolyDataReader, [880](#)
- RequestInformation\_HemodynamicWaveformStorage
  - vtkGDCMPolyDataReader, [880](#)
- RequestInformation\_RTStructureSetStorage
  - vtkGDCMPolyDataReader, [880](#)
- RequestInformationCompat
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [868](#)
- RequestPaddedCompositePixelCode
  - gdcm::ImageCodec, [435](#)
- RequestPlanarConfiguration
  - gdcm::ImageCodec, [435](#)
- Rescale
  - gdcm::Rescaler, [636](#)
- RescaleFunctionIntoBestFit
  - gdcm::Rescaler, [636](#)
- Rescaler
  - gdcm::Rescaler, [635](#)
- ReserveDataElement
  - gdcm::FileStreamer, [394](#)
- ReserveGroupDataElement
  - gdcm::FileStreamer, [395](#)
- reset
  - gdcm::terminal, [137](#)
- ResetHandledDataSet
  - gdcm::network::ULConnectionCallback, [822](#)
- RetrieveSOPInstanceUIDFromIndex
  - gdcm::DirectoryHelper, [332](#)
- RetrieveSOPInstanceUIDFromZPosition
  - gdcm::DirectoryHelper, [332](#)
- reverse
  - gdcm::terminal, [137](#)
- reverse\_iterator
  - gdcm::CodeString, [247](#)
  - gdcm::LO, [495](#)
  - gdcm::String, [709](#)
- RoleSelectionSub
  - gdcm::network::RoleSelectionSub, [641](#)
- SAGITTAL
  - gdcm::Orientation, [545](#)
- SH
  - gdcm::VR, [853](#)
- SHA1
  - gdcm::SHA1, [680](#)
- SHComp
  - gdcm, [124](#)
- SIEMENS
  - gdcm::Dicts, [323](#)
- SINGLEBIT
  - gdcm::PixelFormat, [574](#)
- SL
  - gdcm::VR, [853](#)
- SLICE\_ORIENTATION\_XY
  - vtkImageColorViewer, [895](#)

- SLICE\_ORIENTATION\_XZ
  - vtkImageColorViewer, [895](#)
- SLICE\_ORIENTATION\_YZ
  - vtkImageColorViewer, [895](#)
- SOPClassExtendedNegociationSub
  - gdcm::network::SOPClassExtendedNegociationSub, [689](#)
- SOPInstanceUID
  - vtkRTStructSetProperties, [916](#)
- SPM2AVG152PDFFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2AVG152T1FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2AVG152T2FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2AVG305T1FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2BRAINMASKFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2CSFFFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2EPIFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2FILT1FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2GRAYFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2PDFFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2PETFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2SINGLESUBJT1FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2SPECTFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2T1FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2T2FrameofReference
  - gdcm::UIDs, [773](#)
- SPM2TRANSMFrameofReference
  - gdcm::UIDs, [773](#)
- SPM2WHITEFrameofReference
  - gdcm::UIDs, [773](#)
- SQ
  - gdcm::VR, [853](#)
- SS
  - gdcm::VR, [853](#)
- ST
  - gdcm::VR, [853](#)
- STATES
  - gdcm::Surface, [719](#)
- STATES\_END
  - gdcm::Surface, [719](#)
- STComp
  - gdcm, [124](#)
- SURFACE
  - gdcm::Surface, [719](#)
- SV10
  - gdcm::CSAHeader, [270](#)
- ScalarType
  - gdcm::PixelFormat, [573](#)
- Scale
  - vtkGDCMImageReader, [865](#)
  - vtkGDCMImageReader2, [871](#)
- Scan
  - gdcm::Scanner, [647](#)
- Scanner
  - gdcm::Scanner, [645](#)
- SecondaryCaptureImageStorage
  - gdcm::MediaStorage, [508](#)
  - gdcm::UIDs, [775](#)
- Segment
  - gdcm::Segment, [650](#)
- SegmentAlgorithmName
  - gdcm::Segment, [651](#)
- SegmentAlgorithmType
  - gdcm::Segment, [651](#)
- SegmentDescription
  - gdcm::Segment, [651](#)
- SegmentLabel
  - gdcm::Segment, [651](#)
- SegmentMap
  - gdcm::SegmentReader, [656](#)
- SegmentNumber
  - gdcm::Segment, [652](#)
- SegmentReader
  - gdcm::SegmentReader, [656](#)
- SegmentVector
  - gdcm::SegmentReader, [656](#)
  - gdcm::SegmentWriter, [658](#)
- SegmentWriter
  - gdcm::SegmentWriter, [658](#)
- Segmentation
  - gdcm::MediaStorage, [510](#)
- SegmentationStorage
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [776](#)
- SegmentedPaletteColorLookupTable
  - gdcm::SegmentedPaletteColorLookupTable, [653](#)
- Segments
  - gdcm::SegmentReader, [656](#)
  - gdcm::SegmentWriter, [659](#)
- Selection
  - gdcm::Sorter, [694](#)
- SelectionMap
  - gdcm::Sorter, [692](#)
- Self
  - gdcm::AnonymizeEvent, [153](#)

- gdcm::DataEvent, [293](#)
- gdcm::DataSetEvent, [303](#)
- gdcm::FileNameEvent, [388](#)
- gdcm::MemberCommand, [514](#)
- gdcm::ProgressEvent, [609](#)
- gdcm::SimpleMemberCommand, [682](#)
- SendEcho
  - gdcm::ServiceClassUser, [677](#)
  - gdcm::network::ULConnectionManager, [826](#)
- SendFind
  - gdcm::ServiceClassUser, [677](#)
  - gdcm::network::ULConnectionManager, [826](#)
- SendMove
  - gdcm::ServiceClassUser, [677](#)
  - gdcm::network::ULConnectionManager, [826](#)
- SendStore
  - gdcm::ServiceClassUser, [677](#), [678](#)
  - gdcm::network::ULConnectionManager, [826](#)
- Separator
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::PersonName, [567](#)
- SequenceLengthField
  - gdcm::SequenceOfItems, [669](#)
- SequenceOfFragments
  - gdcm::SequenceOfFragments, [661](#)
- SequenceOfItems
  - gdcm::SequenceOfItems, [667](#)
- SerieHelper
  - gdcm::SerieHelper, [671](#)
- SerieRestrictions
  - gdcm::SerieHelper, [671](#)
- Series
  - gdcm::Series, [673](#)
- SeriesInstanceUID
  - vtkRTStructSetProperties, [916](#)
- ServiceClassApplicationInformation
  - gdcm::network::ServiceClassApplicationInformation, [673](#)
- ServiceClassUser
  - gdcm::ServiceClassUser, [676](#)
- Set
  - gdcm::Attribute, [174](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [179](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [186](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
- SetAETitle
  - gdcm::ServiceClassUser, [678](#)
- SetAbstractSyntax
  - gdcm::PresentationContext, [594](#)
  - gdcm::network::PresentationContextRQ, [599](#)
- SetAlgorithmFamily
  - gdcm::Surface, [721](#)
- SetAlgorithmName
  - gdcm::Surface, [721](#)
- SetAlgorithmVersion
  - gdcm::Surface, [721](#)
- SetAnatomicRegion
  - gdcm::Segment, [651](#)
- SetArray
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
- SetAxisOfRotation
  - gdcm::Surface, [721](#)
- SetBitPosition
  - gdcm::Overlay, [551](#)
- SetBitSample
  - gdcm::JPEGCodec, [486](#)
- SetBitsAllocated
  - gdcm::Overlay, [551](#)
  - gdcm::PixelFormat, [576](#)
- SetBitsStored
  - gdcm::PixelFormat, [576](#)
- SetBlob
  - gdcm::ApplicationEntity, [162](#)
  - gdcm::PersonName, [566](#)
  - gdcm::network::PresentationDataValue, [600](#)
- SetBlueLUT
  - gdcm::LookupTable, [500](#)
- SetBufferLength
  - gdcm::JPEGLSCodec, [490](#)
  - gdcm::PNMCodec, [590](#)
  - gdcm::RLECodec, [640](#)
- SetByteSwapTag
  - gdcm::ByteSwapFilter, [224](#)
- SetByteValue
  - gdcm::Attribute, [174](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [179](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [186](#)
  - gdcm::CSAElement, [267](#)
  - gdcm::DataElement, [289](#)
- SetByteValueNoSwap
  - gdcm::Attribute, [174](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [179](#)
- SetCallbackFunction
  - gdcm::MemberCommand, [515](#)
  - gdcm::SimpleMemberCommand, [683](#)
- SetCalledAETitle
  - gdcm::ServiceClassUser, [678](#)
  - gdcm::network::AAssociateACPDU, [143](#)
  - gdcm::network::AAssociateRQPDU, [148](#)
- SetCallingAETitle
  - gdcm::network::AAssociateACPDU, [143](#)
  - gdcm::network::AAssociateRQPDU, [148](#)

- SetCenterOfRotation
  - gdcm::Surface, [721](#)
- SetChangePrivateTags
  - gdcm::FileExplicitFilter, [378](#)
- SetCheckFileMetaInformation
  - gdcm::Writer, [920](#)
- SetCipherType
  - gdcm::CAPICryptographicMessageSyntax, [232](#)
  - gdcm::CryptographicMessageSyntax, [263](#)
  - gdcm::OpenSSLCryptographicMessageSyntax, [540](#)
  - gdcm::OpenSSL7CryptographicMessageSyntax, [544](#)
- SetColor
  - gdcm::Printer, [604](#)
- SetColorLevel
  - vtkImageColorViewer, [896](#)
- SetColorWindow
  - vtkImageColorViewer, [896](#)
- SetColumns
  - gdcm::Bitmap, [214](#)
  - gdcm::Overlay, [551](#)
- SetCommand
  - gdcm::network::PresentationDataValue, [600](#)
- SetComponents
  - gdcm::PersonName, [566](#)
- SetCompressIconImage
  - gdcm::ImageChangeTransferSyntax, [427](#)
- SetComputeZSpacing
  - gdcm::IPPSorter, [466](#)
- SetCoordinateStartValue
  - gdcm::Curve, [282](#)
- SetCoordinateStepValue
  - gdcm::Curve, [282](#)
- SetCryptographicMessageSyntax
  - gdcm::Anonymizer, [158](#)
- SetCurve
  - gdcm::Curve, [282](#)
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageReader2, [869](#)
- SetCurveDataDescriptor
  - gdcm::Curve, [282](#)
- SetCurveDescription
  - gdcm::Curve, [282](#)
- SetData
  - gdcm::DataEvent, [294](#)
- SetDataElement
  - gdcm::Bitmap, [214](#)
- SetDataSet
  - gdcm::File, [368](#)
  - gdcm::network::PresentationDataValue, [600](#)
- SetDataSetTransferSyntax
  - gdcm::FileMetaInformation, [383](#)
- SetDataValueRepresentation
  - gdcm::Curve, [282](#)
- SetDebug
  - gdcm::Trace, [756](#)
- SetDebugStream
  - gdcm::Trace, [756](#)
- SetDefaultTransferSyntax
  - gdcm::PresentationContextGenerator, [597](#)
- SetDerivationCodeSequenceCodeValue
  - gdcm::FileDerivation, [376](#)
- SetDerivationDescription
  - gdcm::FileDerivation, [376](#)
- SetDescription
  - gdcm::CSAHeaderDictEntry, [275](#)
  - gdcm::ModuleEntry, [525](#)
  - gdcm::Overlay, [551](#)
- SetDescriptor
  - gdcm::DICOMDIRGenerator, [313](#)
- SetDictName
  - gdcm::DictConverter, [317](#)
- SetDicts
  - gdcm::PythonFilter, [613](#)
  - gdcm::StringFilter, [712](#)
- SetDimension
  - gdcm::Bitmap, [214](#)
- SetDimensions
  - gdcm::Bitmap, [215](#)
  - gdcm::Curve, [282](#)
  - gdcm::ImageCodec, [433](#)
- SetDimensionsValue
  - gdcm::ImageHelper, [442](#)
- SetDirectionCosines
  - gdcm::Image, [416](#)
  - vtkGDCMImageWriter, [874](#)
- SetDirectionCosinesFromImageOrientationPatient
  - vtkGDCMImageWriter, [874](#)
- SetDirectionCosinesTolerance
  - gdcm::IPPSorter, [466](#)
- SetDirectionCosinesValue
  - gdcm::ImageHelper, [442](#)
- SetDirectory
  - gdcm::SerieHelper, [672](#)
  - gdcm::network::ULWritingCallback, [829](#)
- SetDisplayId
  - vtkImageColorViewer, [896](#)
- SetDomain
  - gdcm::BoxRegion, [221](#)
- SetDropDuplicatePositions
  - gdcm::IPPSorter, [466](#)
- SetElement
  - gdcm::Tag, [747](#)
- SetElementHandler
  - gdcm::Parser, [556](#)
- SetElementTag
  - gdcm::Tag, [747](#)
- SetElementXX

- gdcm::DictEntry, [320](#)
- SetError
  - gdcm::Trace, [756](#)
- SetErrorStream
  - gdcm::Trace, [756](#)
- SetEvent
  - gdcm::network::ULEvent, [827](#)
- SetFile
  - gdcm::Anonymizer, [158](#)
  - gdcm::DICOMDIRGenerator, [313](#)
  - gdcm::FileDerivation, [376](#)
  - gdcm::FileExplicitFilter, [378](#)
  - gdcm::IconImageFilter, [410](#)
  - gdcm::Printer, [604](#)
  - gdcm::PythonFilter, [613](#)
  - gdcm::Reader, [631](#)
  - gdcm::SplitMosaicFilter, [697](#)
  - gdcm::StreamImageWriter, [705](#)
  - gdcm::StringFilter, [712](#)
  - gdcm::Validate, [840](#)
  - gdcm::Writer, [920](#)
  - gdcm::XMLPrinter, [925](#)
- SetFileName
  - gdcm::FileNameEvent, [388](#)
  - gdcm::Reader, [631](#)
  - gdcm::StreamImageReader, [701](#)
  - gdcm::StreamImageWriter, [705](#)
  - gdcm::Writer, [920](#)
  - vtkGDCMThreadedImageReader2, [890](#)
- SetFileNames
  - vtkGDCMImageReader, [862](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMThreadedImageReader2, [890](#)
- SetFilePattern
  - vtkGDCMImageReader, [863](#)
  - vtkGDCMImageReader2, [869](#)
- SetFilePrefix
  - vtkGDCMImageReader, [863](#)
  - vtkGDCMImageReader2, [869](#)
- SetFilename
  - gdcm::TableReader, [740](#)
- SetFilenames
  - gdcm::DICOMDIRGenerator, [313](#)
- SetFiles
  - gdcm::FileSet, [392](#)
- SetFiniteVolume
  - gdcm::Surface, [721](#)
- SetForce
  - gdcm::ImageChangeTransferSyntax, [428](#)
  - gdcm::ImageFragmentSplitter, [439](#)
- SetForcePixelSpacing
  - gdcm::ImageHelper, [442](#)
- SetForceRescaleInterceptSlope
  - gdcm::ImageHelper, [442](#)
- SetFragmentSizeMax
  - gdcm::ImageFragmentSplitter, [439](#)
- SetFrameOrigin
  - gdcm::Overlay, [551](#)
- SetFromDataElement
  - gdcm::Attribute, [174](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [179](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [186](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
- SetFromDataSet
  - gdcm::Attribute, [175](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1 >, [180](#)
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, [186](#)
  - gdcm::MediaStorage, [512](#)
- SetFromFile
  - gdcm::MediaStorage, [512](#)
- SetFromHeader
  - gdcm::MediaStorage, [512](#)
- SetFromModality
  - gdcm::MediaStorage, [512](#)
- SetFromSourceImageSequence
  - gdcm::MediaStorage, [512](#)
- SetFromString
  - gdcm::DirectionCosines, [328](#)
- SetFromUID
  - gdcm::UIDs, [786](#)
- SetGreenLUT
  - gdcm::LookupTable, [500](#)
- SetGroup
  - gdcm::Curve, [282](#)
  - gdcm::Overlay, [551](#)
  - gdcm::Tag, [747](#)
- SetGroupXX
  - gdcm::DictEntry, [320](#)
- SetHeader
  - gdcm::File, [368](#)
- SetHighBit
  - gdcm::PixelFormat, [576](#)
- SetHostname
  - gdcm::ServiceClassUser, [678](#)
- SetIE
  - gdcm::IODEntry, [462](#)
- SetIconImage
  - gdcm::Pixmap, [579](#)
- SetImage
  - gdcm::PixmapWriter, [587](#)
  - gdcm::SplitMosaicFilter, [697](#)
- SetImplementationClassUID
  - gdcm::FileMetaInformation, [383](#)

- SetImplementationVersionName
  - gdcm::FileMetaInformation, [383](#)
- SetImplicitFlag
  - gdcm::network::ULConnectionCallback, [822](#)
- SetInput
  - gdcm::BitmapToBitmapFilter, [218](#)
  - gdcm::ImageConverter, [436](#)
  - vtkImageColorViewer, [896](#)
- SetInputConnection
  - vtkImageColorViewer, [896](#)
- SetInputFileName
  - gdcm::DictConverter, [317](#)
  - gdcm::FileAnonymizer, [371](#)
  - gdcm::FileChangeTransferSyntax, [374](#)
- SetIntercept
  - gdcm::Image, [416](#)
  - gdcm::Rescaler, [636](#)
- SetKey
  - gdcm::CSAElement, [267](#)
- SetKeyword
  - gdcm::DictEntry, [320](#)
- SetLUT
  - gdcm::Bitmap, [215](#)
  - gdcm::ImageCodec, [434](#)
  - gdcm::LookupTable, [500](#)
  - gdcm::SegmentedPaletteColorLookupTable, [653](#)
- SetLastElement
  - gdcm::ParseException, [553](#)
- SetLastFragment
  - gdcm::network::PresentationDataValue, [601](#)
- SetLength
  - gdcm::ByteValue, [229](#)
  - gdcm::Element< TVR, VM::VM1\_2 >, [339](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
  - gdcm::Element< TVR, VM::VM2\_2n >, [343](#)
  - gdcm::Element< TVR, VM::VM2\_n >, [344](#)
  - gdcm::Element< TVR, VM::VM3\_3n >, [346](#)
  - gdcm::Element< TVR, VM::VM3\_n >, [347](#)
  - gdcm::RLECodec, [640](#)
  - gdcm::SequenceOfFragments, [663](#)
  - gdcm::SequenceOfItems, [669](#)
  - gdcm::Value, [842](#)
- SetLengthOnly
  - gdcm::ByteValue, [229](#)
  - gdcm::Value, [842](#)
- SetLengthToUndefined
  - gdcm::SequenceOfItems, [669](#)
- SetLoadMode
  - gdcm::SerieHelper, [672](#)
- SetLookupTable
  - vtkImageMapToColors16, [901](#)
- SetLossless
  - gdcm::JPEGCodec, [486](#)
  - gdcm::JPEGLSCodec, [490](#)
- SetLossyError
  - gdcm::JPEGLSCodec, [490](#)
- SetLossyFlag
  - gdcm::Bitmap, [215](#)
  - gdcm::ImageCodec, [434](#)
  - gdcm::PVRGCodec, [612](#)
- SetManifold
  - gdcm::Surface, [721](#)
- SetMaxPDULength
  - gdcm::network::ULConnectionInfo, [823](#)
- SetMaxPDUSize
  - gdcm::network::ULConnection, [820](#)
- SetMaximumLength
  - gdcm::network::MaximumLengthSub, [504](#)
- SetMaximumPointDistance
  - gdcm::Surface, [721](#)
- SetMeanPointDistance
  - gdcm::Surface, [721](#)
- SetMedicalImageProperties
  - vtkGDCMImageReader, [863](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMPolyDataWriter, [882](#)
- SetMergeModeToAbstractSyntax
  - gdcm::PresentationContextGenerator, [597](#)
- SetMergeModeToTransferSyntax
  - gdcm::PresentationContextGenerator, [597](#)
- SetMeshPrimitive
  - gdcm::Surface, [721](#)
- SetMessageHeader
  - gdcm::network::PresentationDataValue, [601](#)
- SetMinMaxForPixelType
  - gdcm::Rescaler, [636](#)
- SetName
  - gdcm::CSAElement, [267](#)
  - gdcm::CSAHeaderDictEntry, [276](#)
  - gdcm::DictEntry, [320](#)
  - gdcm::IODEntry, [462](#)
  - gdcm::Macro, [502](#)
  - gdcm::Module, [522](#)
  - gdcm::ModuleEntry, [525](#)
  - gdcm::PDBelement, [560](#)
  - gdcm::network::AbstractSyntax, [151](#)
  - gdcm::network::ApplicationContext, [161](#)
  - gdcm::network::TransferSyntaxSub, [762](#)
- SetNameFromUID
  - gdcm::network::AbstractSyntax, [151](#)
  - gdcm::network::TransferSyntaxSub, [762](#)
- SetNeedByteSwap
  - gdcm::Bitmap, [215](#)
  - gdcm::ImageCodec, [434](#)
- SetNeedOverlayCleanup
  - gdcm::ImageCodec, [434](#)
- SetNestedDataSet



- gdcm::Item, 470
- SetNoOfItems
  - gdcm::CSAElement, 267
- SetNoSwap
  - gdcm::Element, 337
  - gdcm::Element< TVR, VM::VM1\_n >, 341
- SetNumberOfCurves
  - gdcm::Pixmap, 579
- SetNumberOfDimensions
  - gdcm::Bitmap, 215
  - gdcm::ImageCodec, 434
- SetNumberOfFilenames
  - gdcm::FilenameGenerator, 390
- SetNumberOfFrames
  - gdcm::Overlay, 551
- SetNumberOfInputPorts
  - vtkGDCMPolyDataWriter, 882
- SetNumberOfItems
  - gdcm::SequenceOfItems, 669
- SetNumberOfOverlays
  - gdcm::Pixmap, 579
- SetNumberOfPoints
  - gdcm::Curve, 282
- SetNumberOfResolutions
  - gdcm::JPEG2000Codec, 479
- SetNumberOfSegments
  - gdcm::SegmentWriter, 658
- SetNumberOfSurfacePoints
  - gdcm::Surface, 721
- SetNumberOfSurfaces
  - gdcm::SurfaceWriter, 729
- SetNumberOfTableValues
  - vtkLookupTable16, 911
- SetNumberOfValues
  - gdcm::Attribute< Group, Element, TVR, VM::VM1\_n >, 186
- SetNumberOfVectors
  - gdcm::Surface, 721
- SetObliquityThresholdCosineValue
  - gdcm::Orientation, 546
- SetOffScreenRendering
  - vtkImageColorViewer, 897
- SetOrigin
  - gdcm::Image, 416
  - gdcm::Overlay, 551
- SetOriginValue
  - gdcm::ImageHelper, 442
- SetOutputDimensions
  - gdcm::IconImageGenerator, 412
- SetOutputFileName
  - gdcm::DictConverter, 317
  - gdcm::FileAnonymizer, 371
  - gdcm::FileChangeTransferSyntax, 374
  - gdcm::FileStreamer, 395
  - SetOutputFormatToLuminance
    - vtkImageMapToColors16, 901
  - SetOutputFormatToLuminanceAlpha
    - vtkImageMapToColors16, 901
  - SetOutputFormatToRGB
    - vtkImageMapToColors16, 901
  - SetOutputFormatToRGBA
    - vtkImageMapToColors16, 901
  - SetOutputType
    - gdcm::DictConverter, 317
  - SetOutsideValuePixel
    - gdcm::IconImageGenerator, 412
  - SetOverlay
    - gdcm::Overlay, 551
  - SetOverlayVisibility
    - vtkImageColorViewer, 897
  - SetOwner
    - gdcm::PrivateTag, 607
  - SetPDU
    - gdcm::network::ULEvent, 827
  - SetParentId
    - vtkImageColorViewer, 897
  - SetPassword
    - gdcm::CAPICryptographicMessageSyntax, 232
    - gdcm::CryptographicMessageSyntax, 263
    - gdcm::OpenSSLCryptographicMessageSyntax, 540
    - gdcm::OpenSSL7CryptographicMessageSyntax, 544
  - SetPattern
    - gdcm::FilenameGenerator, 390
  - SetPermissions
    - gdcm::System, 736
  - SetPhotometricInterpretation
    - gdcm::Bitmap, 215
    - gdcm::ImageChangePhotometricInterpretation, 421
    - gdcm::ImageCodec, 434
  - SetPixelFormat
    - gdcm::Bitmap, 215
    - gdcm::ImageCodec, 434
    - gdcm::JPEGCodec, 486
    - gdcm::Rescaler, 636
  - SetPixelMinMax
    - gdcm::IconImageGenerator, 412
  - SetPixelRepresentation
    - gdcm::PixelFormat, 576
  - SetPixmap
    - gdcm::IconImageGenerator, 412
    - gdcm::PixmapWriter, 588
  - SetPlanarConfiguration
    - gdcm::Bitmap, 215
    - gdcm::ImageChangePlanarConfiguration, 424
    - gdcm::ImageCodec, 434
  - SetPointCoordinatesData
    - gdcm::Surface, 721

- SetPointPositionAccuracy
  - gdcm::Surface, [722](#)
- SetPointsBoundingBoxCoordinates
  - gdcm::Surface, [722](#)
- SetPort
  - gdcm::ServiceClassUser, [678](#)
- SetPortSCP
  - gdcm::ServiceClassUser, [678](#)
- SetPosition
  - vtkImageColorViewer, [897](#)
- SetPreamble
  - gdcm::FileMetaInformation, [383](#)
- SetPrefix
  - gdcm::FilenameGenerator, [390](#)
- SetPresentationContextID
  - gdcm::PresentationContext, [594](#)
  - gdcm::network::PresentationContextAC, [595](#)
  - gdcm::network::PresentationContextRQ, [599](#)
  - gdcm::network::PresentationDataValue, [601](#)
- SetPresentationContexts
  - gdcm::ServiceClassUser, [678](#)
  - gdcm::network::ULConnection, [820](#)
- SetPrettyPrint
  - gdcm::JSON, [491](#)
- SetPrimitiveData
  - gdcm::MeshPrimitive, [519](#)
- SetPrimitiveType
  - gdcm::MeshPrimitive, [519](#)
- SetPrimitivesData
  - gdcm::MeshPrimitive, [519](#)
- SetPrivateCreator
  - gdcm::Tag, [747](#)
- SetProcessingAlgorithm
  - gdcm::Surface, [722](#)
- SetProgress
  - gdcm::ProgressEvent, [609](#)
- SetPropertyCategory
  - gdcm::Segment, [651](#)
- SetPropertyType
  - gdcm::Segment, [651](#)
- SetPurposeOfReferenceCodeSequenceCodeValue
  - gdcm::FileDerivation, [376](#)
- SetQuality
  - gdcm::JPEG2000Codec, [479](#)
  - gdcm::JPEGCodec, [486](#)
- SetRTStructSetProperties
  - vtkGDCMPolyDataWriter, [883](#)
- SetRate
  - gdcm::JPEG2000Codec, [479](#)
- SetReason
  - gdcm::network::AAAbortPDU, [140](#)
  - gdcm::network::PresentationContextAC, [595](#)
- SetRecommendedDisplayCIELabValue
  - gdcm::Surface, [722](#)
- SetRecommendedDisplayGrayscaleValue
  - gdcm::Surface, [722](#)
- SetRecommendedPresentationOpacity
  - gdcm::Surface, [722](#)
- SetRecommendedPresentationType
  - gdcm::Surface, [722](#)
- SetRecomputeItemLength
  - gdcm::FileExplicitFilter, [378](#)
- SetRecomputeSequenceLength
  - gdcm::FileExplicitFilter, [378](#)
- SetRedLUT
  - gdcm::LookupTable, [500](#)
- SetRef
  - gdcm::IODEntry, [462](#)
- SetRegion
  - gdcm::ImageRegionReader, [449](#)
- SetRenderWindow
  - vtkImageColorViewer, [897](#)
- SetRenderer
  - vtkImageColorViewer, [897](#)
- SetRescaleInterceptSlopeValue
  - gdcm::ImageHelper, [442](#)
- SetRetired
  - gdcm::DictEntry, [320](#)
- SetReversible
  - gdcm::JPEG2000Codec, [479](#)
- SetRoot
  - gdcm::UIDGenerator, [767](#)
- SetRootDirectory
  - gdcm::DICOMDIRGenerator, [313](#)
- SetRows
  - gdcm::Bitmap, [215](#)
  - gdcm::Overlay, [552](#)
- SetSamplesPerPixel
  - gdcm::PixelFormat, [576](#)
- SetScalarType
  - gdcm::PixelFormat, [576](#)
- SetSearchParameter
  - gdcm::BaseRootQuery, [203](#)
- SetSegmentAlgorithmName
  - gdcm::Segment, [651](#)
- SetSegmentAlgorithmType
  - gdcm::Segment, [651](#)
- SetSegmentDescription
  - gdcm::Segment, [651](#)
- SetSegmentLabel
  - gdcm::Segment, [651](#)
- SetSegmentNumber
  - gdcm::Segment, [651](#)
- SetSegments
  - gdcm::SegmentWriter, [658](#)
- SetSize
  - vtkImageColorViewer, [897](#)
- SetSlice



- vtkImageColorViewer, [897](#)
- SetSliceOrientation
  - vtkImageColorViewer, [897](#)
- SetSliceOrientationToXY
  - vtkImageColorViewer, [897](#)
- SetSliceOrientationToXZ
  - vtkImageColorViewer, [897](#)
- SetSliceOrientationToYZ
  - vtkImageColorViewer, [897](#)
- SetSlope
  - gdcm::Image, [416](#)
  - gdcm::Rescaler, [636](#)
- SetSortFunction
  - gdcm::Sorter, [693](#)
- SetSource
  - gdcm::network::AAAbortPDU, [141](#)
- SetSourceApplicationEntityTitle
  - gdcm::FileMetaInformation, [383](#)
- SetSpacing
  - gdcm::Image, [416](#)
- SetSpacingValue
  - gdcm::ImageHelper, [442](#)
- SetState
  - gdcm::network::ULConnection, [820](#)
- SetStream
  - gdcm::Reader, [632](#)
  - gdcm::StreamImageReader, [702](#)
  - gdcm::StreamImageWriter, [705](#)
  - gdcm::Trace, [756](#)
  - gdcm::Writer, [921](#)
- SetStreamToFile
  - gdcm::Trace, [756](#)
- SetStyle
  - gdcm::Printer, [604](#)
  - gdcm::XMLPrinter, [925](#)
- SetSurfaceComments
  - gdcm::Surface, [722](#)
- SetSurfaceCount
  - gdcm::Segment, [651](#)
- SetSurfaceNumber
  - gdcm::Surface, [722](#)
- SetSurfaceProcessing
  - gdcm::Surface, [722](#)
- SetSurfaceProcessingDescription
  - gdcm::Surface, [722](#)
- SetSurfaceProcessingRatio
  - gdcm::Surface, [722](#)
- SetSyngoDT
  - gdcm::CSAElement, [267](#)
- SetTag
  - gdcm::AnonymizeEvent, [153](#)
  - gdcm::DataElement, [289](#)
- SetTargetPixelType
  - gdcm::Rescaler, [636](#)
- SetTemplateFileName
  - gdcm::FileStreamer, [395](#)
- SetTileSize
  - gdcm::JPEG2000Codec, [479](#)
- SetTimeout
  - gdcm::ServiceClassUser, [678](#)
  - gdcm::network::ARTIMTimer, [167](#)
- SetToUndefined
  - gdcm::VL, [846](#)
- SetTransferSyntax
  - gdcm::Bitmap, [216](#)
  - gdcm::FileChangeTransferSyntax, [374](#)
  - gdcm::ImageChangeTransferSyntax, [428](#)
  - gdcm::network::PresentationContextAC, [595](#)
- SetTuple
  - gdcm::network::RoleSelectionSub, [641](#)
  - gdcm::network::SOPClassExtendedNegociationSub, [689](#)
  - gdcm::network::ServiceClassApplicationInformation, [673](#)
- SetType
  - gdcm::ModuleEntry, [525](#)
  - gdcm::Overlay, [552](#)
- SetTypeOfData
  - gdcm::Curve, [282](#)
- SetUsage
  - gdcm::IODEntry, [462](#)
- SetUseSeriesDetails
  - gdcm::SerieHelper, [672](#)
- SetUseTargetPixelType
  - gdcm::Rescaler, [636](#)
- SetUseVRUN
  - gdcm::FileExplicitFilter, [379](#)
- SetUserCodec
  - gdcm::ImageChangeTransferSyntax, [428](#)
- SetUserData
  - gdcm::Parser, [556](#)
- SetUserInformation
  - gdcm::network::AAAssociateRQPDU, [148](#)
- SetVL
  - gdcm::DataElement, [290](#)
- SetVLToUndefined
  - gdcm::DataElement, [290](#)
- SetVM
  - gdcm::CSAElement, [267](#)
  - gdcm::CSAHeaderDictEntry, [276](#)
  - gdcm::DictEntry, [320](#)
- SetVR
  - gdcm::CSAElement, [267](#)
  - gdcm::CSAHeaderDictEntry, [276](#)
  - gdcm::DataElement, [290](#)
  - gdcm::DictEntry, [320](#)
- SetValue
  - gdcm::Attribute, [175](#)

- gdcmm::Attribute< Group, Element, TVR, VM::VM1 >, 180
- gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, 186
- gdcmm::CSAElement, 267
- gdcmm::DataElement, 289
- gdcmm::Element, 337
- gdcmm::Element< TVR, VM::VM1\_n >, 341
- gdcmm::PDSElement, 560
- SetValueFieldLength
  - gdcmm::DataElement, 290
- SetValues
  - gdcmm::Attribute, 175
  - gdcmm::Attribute< Group, Element, TVR, VM::VM1\_n >, 187
- SetVectorAccuracy
  - gdcmm::Surface, 722
- SetVectorCoordinateData
  - gdcmm::Surface, 722
- SetVectorDimensionality
  - gdcmm::Surface, 722
- SetWarning
  - gdcmm::Trace, 756
- SetWarningStream
  - gdcmm::Trace, 756
- SetWindowId
  - vtkImageColorViewer, 898
- SetWriteDataSetOnly
  - gdcmm::Writer, 921
- SetZSpacingTolerance
  - gdcmm::IPPSorter, 466
- setattribute
  - gdcmm::terminal, 138
- setbgcolor
  - gdcmm::terminal, 138
- setfgcolor
  - gdcmm::terminal, 138
- setmode
  - gdcmm::terminal, 138
- SetupInteractor
  - vtkImageColorViewer, 897
- Shift
  - vtkGDCMImageReader, 865
  - vtkGDCMImageReader2, 871
- ShiftEnd
  - gdcmm::ByteBuffer, 222
- ShowAbort
  - gdcmm::SimpleSubjectWatcher, 685
- ShowAnonymization
  - gdcmm::SimpleSubjectWatcher, 685
- ShowData
  - gdcmm::SimpleSubjectWatcher, 685
- ShowDataSet
  - gdcmm::SimpleSubjectWatcher, 685
- ShowFileName
  - gdcmm::SimpleSubjectWatcher, 685
- ShowIteration
  - gdcmm::SimpleSubjectWatcher, 685
- ShowProgress
  - gdcmm::SimpleSubjectWatcher, 685
- SimpleMemberCommand
  - gdcmm::SimpleMemberCommand, 683
- SimpleSubjectWatcher
  - gdcmm::SimpleSubjectWatcher, 684
- SingleSerieUIDFileSetHT
  - gdcmm::SerieHelper, 672
- SingleSerieUIDFileSetmap
  - gdcmm::SerieHelper, 671
- Size
  - gdcmm::CodeString, 248
  - gdcmm::DataSet, 301
  - gdcmm::GroupDict, 408
  - gdcmm::network::AAAbortPDU, 141
  - gdcmm::network::AAAssociateACPDU, 143
  - gdcmm::network::AAAssociateRJPDU, 145
  - gdcmm::network::AAAssociateRQPDU, 149
  - gdcmm::network::AResetRPPDU, 164
  - gdcmm::network::AResetRQPDU, 166
  - gdcmm::network::AbstractSyntax, 151
  - gdcmm::network::ApplicationContext, 161
  - gdcmm::network::AsynchronousOperationsWindow←Sub, 169
  - gdcmm::network::BasePDU, 200
  - gdcmm::network::ImplementationClassUIDSub, 454
  - gdcmm::network::ImplementationVersionNameSub, 455
  - gdcmm::network::MaximumLengthSub, 504
  - gdcmm::network::PDataTFPDU, 558
  - gdcmm::network::PresentationContextAC, 595
  - gdcmm::network::PresentationContextRQ, 599
  - gdcmm::network::PresentationDataValue, 601
  - gdcmm::network::RoleSelectionSub, 641
  - gdcmm::network::SOPClassExtendedNegotiationSub, 689
  - gdcmm::network::ServiceClassApplicationInformation, 674
  - gdcmm::network::TransferSyntaxSub, 762
  - gdcmm::network::UserInformation, 838
- size\_type
  - gdcmm::CodeString, 247
  - gdcmm::LO, 495
  - gdcmm::String, 709
- SizeType
  - gdcmm::DataSet, 297
  - gdcmm::FilenameGenerator, 389
  - gdcmm::IOD, 459
  - gdcmm::NestedModuleEntries, 533
  - gdcmm::PresentationContext, 593

- gdcm::PresentationContextGenerator, 596
- gdcm::SequenceOfFragments, 661
- gdcm::SequenceOfItems, 667
- gdcm::network::AAAssociateACPDU, 143
- gdcm::network::AAAssociateRQPDU, 147
- gdcm::network::PDataTFPDU, 558
- gdcm::network::PresentationContextRQ, 598
- Slice
  - vtkImageColorViewer, 899
- SliceOrientation
  - vtkImageColorViewer, 899
- SmartPointer
  - gdcm::Object, 536
  - gdcm::SmartPointer, 687
- Sort
  - gdcm::IPPSorter, 466
  - gdcm::Sorter, 693
- SortFunc
  - gdcm::Sorter, 694
- SortFunction
  - gdcm::Sorter, 692
- Sorter
  - gdcm::Sorter, 693
- SpacialFiducialsStorage
  - gdcm::MediaStorage, 509
- SpacialRegistrationStorage
  - gdcm::MediaStorage, 509
- Spacing
  - gdcm::Spacing, 695
- SpacingType
  - gdcm::Spacing, 695
- SpatialFiducialsStorage
  - gdcm::UIDs, 776
- SpatialRegistrationStorage
  - gdcm::UIDs, 776
- Spectroscopy
  - gdcm::Spectroscopy, 696
- Split
  - gdcm::ImageFragmentSplitter, 439
  - gdcm::SplitMosaicFilter, 697
- SplitExtent
  - vtkGDCMThreadedImageReader2, 890
- SplitMosaicFilter
  - gdcm::SplitMosaicFilter, 697
- Squeeze
  - gdcm::ApplicationEntity, 162
- StableSort
  - gdcm::Sorter, 693
- StandaloneCurveStorage
  - gdcm::MediaStorage, 509
- StandaloneCurveStorageRetired
  - gdcm::UIDs, 775
- StandaloneModalityLUTStorage
  - gdcm::MediaStorage, 509
- StandaloneModalityLUTStorageRetired
  - gdcm::UIDs, 776
- StandaloneOverlayStorage
  - gdcm::MediaStorage, 509
- StandaloneOverlayStorageRetired
  - gdcm::UIDs, 775
- StandalonePETCurveStorageRetired
  - gdcm::UIDs, 777
- StandaloneVOILUTStorage
  - gdcm::MediaStorage, 509
- StandaloneVOILUTStorageRetired
  - gdcm::UIDs, 776
- Start
  - gdcm::network::ARTIMTimer, 167
- StartAssociation
  - gdcm::ServiceClassUser, 679
- StartDataElement
  - gdcm::FileStreamer, 395
- StartElement
  - gdcm::TableReader, 740
  - gdcm::XMLDictReader, 923
  - gdcm::XMLPrivateDictReader, 927
- StartElementHandler
  - gdcm::Parser, 555
- StartEncode
  - gdcm::ImageCodec, 434
  - gdcm::JPEG2000Codec, 479
  - gdcm::JPEGCodec, 486
  - gdcm::JPEGLSCodec, 490
  - gdcm::RLECodec, 640
- StartFilter
  - gdcm::SimpleSubjectWatcher, 685
- StartGroupDataElement
  - gdcm::FileStreamer, 395
- StereometricRelationshipStorage
  - gdcm::UIDs, 776
- Stop
  - gdcm::network::ARTIMTimer, 167
- StopAssociation
  - gdcm::ServiceClassUser, 679
- StopDataElement
  - gdcm::FileStreamer, 395
- StopEncode
  - gdcm::ImageCodec, 434
  - gdcm::JPEG2000Codec, 479
  - gdcm::JPEGCodec, 486
  - gdcm::JPEGLSCodec, 490
  - gdcm::RLECodec, 640
- StopGroupDataElement
  - gdcm::FileStreamer, 395
- StopProtocol
  - gdcm::network::ULConnection, 820
- StorageCommitmentPullModelSOPClassRetired
  - gdcm::UIDs, 774

- StorageCommitmentPullModelSOPInstanceRetired
  - gdcm::UIDs, [774](#)
- StorageCommitmentPushModelSOPClass
  - gdcm::UIDs, [774](#)
- StorageCommitmentPushModelSOPInstance
  - gdcm::UIDs, [774](#)
- StorageServiceClass
  - gdcm::UIDs, [774](#)
- StoredPrintStorageSOPClassRetired
  - gdcm::UIDs, [775](#)
- StrCaseCmp
  - gdcm::System, [736](#)
- StrNCaseCmp
  - gdcm::System, [736](#)
- StrSep
  - gdcm::System, [736](#)
- StrTokR
  - gdcm::System, [736](#)
- Stream
  - gdcm::Writer, [921](#)
- StreamImageReader
  - gdcm::Reader, [632](#)
  - gdcm::StreamImageReader, [700](#)
- StreamImageWriter
  - gdcm::StreamImageWriter, [704](#)
  - gdcm::Writer, [921](#)
- String
  - gdcm::String, [709](#), [710](#)
- StringFilter
  - gdcm::StringFilter, [711](#)
- StructureSetDate
  - vtkRTStructSetProperties, [916](#)
- StructureSetLabel
  - vtkRTStructSetProperties, [916](#)
- StructureSetName
  - vtkRTStructSetProperties, [916](#)
- StructureSetTime
  - vtkRTStructSetProperties, [916](#)
- Study
  - gdcm::Study, [713](#)
- StudyComponentManagementSOPClass
  - gdcm::MediaStorage, [509](#)
- StudyComponentManagementSOPClassRetired
  - gdcm::UIDs, [774](#)
- StudyInstanceUID
  - vtkRTStructSetProperties, [916](#)
- StudyRootQueryRetrieveInformationModelFIND
  - gdcm::UIDs, [777](#)
- StudyRootQueryRetrieveInformationModelGET
  - gdcm::UIDs, [777](#)
- StudyRootQueryRetrieveInformationModelMOVE
  - gdcm::UIDs, [777](#)
- Subject
  - gdcm::Subject, [715](#)
- SubstanceAdministrationLoggingSOPClass
  - gdcm::UIDs, [774](#)
- SubstanceAdministrationLoggingSOPInstance
  - gdcm::UIDs, [774](#)
- SubstanceApprovalQuerySOPClass
  - gdcm::UIDs, [778](#)
- Superclass
  - gdcm::AnonymizeEvent, [153](#)
  - gdcm::DataEvent, [293](#)
  - gdcm::DataSetEvent, [303](#)
  - gdcm::FileNameEvent, [388](#)
  - gdcm::LO, [495](#)
  - gdcm::ProgressEvent, [609](#)
- Surface
  - gdcm::Surface, [719](#)
- SurfaceCount
  - gdcm::Segment, [652](#)
- SurfaceReader
  - gdcm::SurfaceReader, [727](#)
- SurfaceSegmentationStorage
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [779](#)
- SurfaceVector
  - gdcm::Segment, [650](#)
- SurfaceWriter
  - gdcm::SurfaceWriter, [729](#)
- Surfaces
  - gdcm::Segment, [652](#)
- Swap
  - gdcm::ByteSwap, [223](#)
  - gdcm::SwapperDoOp, [731](#)
  - gdcm::SwapperNoOp, [732](#)
- SwapArray
  - gdcm::SwapperDoOp, [731](#)
  - gdcm::SwapperNoOp, [732](#)
- SwapCode
  - gdcm::SwapCode, [731](#)
- SwapCodeType
  - gdcm::SwapCode, [730](#)
- SwapFromSwapCodeIntoSystem
  - gdcm::ByteSwap, [223](#)
- SwapRange
  - gdcm::ByteSwap, [223](#)
- SwapRangeFromSwapCodeIntoSystem
  - gdcm::ByteSwap, [223](#)
- SyngoDTField
  - gdcm::CSAElement, [268](#)
- SyntaxError
  - gdcm::Parser, [555](#)
- SystemIsBigEndian
  - gdcm::ByteSwap, [223](#)
- SystemIsLittleEndian
  - gdcm::ByteSwap, [223](#)

- T1
  - gdcm::Type, [765](#)
- T1C
  - gdcm::Type, [765](#)
- T2
  - gdcm::Type, [765](#)
- T2C
  - gdcm::Type, [765](#)
- T3
  - gdcm::Type, [765](#)
- TConstMemberFunctionPointer
  - gdcm::MemberCommand, [514](#)
- TM
  - gdcm::VR, [853](#)
- TMComp
  - gdcm, [124](#)
- TMemberFunctionPointer
  - gdcm::MemberCommand, [514](#)
  - gdcm::SimpleMemberCommand, [682](#)
- TRIANGLE
  - gdcm::MeshPrimitive, [518](#)
- TRIANGLE\_FAN
  - gdcm::MeshPrimitive, [518](#)
- TRIANGLE\_STRIP
  - gdcm::MeshPrimitive, [518](#)
- TS
  - gdcm::Bitmap, [217](#)
- TS\_END
  - gdcm::TransferSyntax, [760](#)
- TSName
  - gdcm::UIDs, [772](#)
- TSType
  - gdcm::TransferSyntax, [759](#)
  - gdcm::UIDs, [779](#)
- TYPETOENCODING
  - gdcm, [130](#)
  - gdcmVR.h, [1187](#)
- TYPETOLENGTH
  - gdcmVM.h, [1185](#)
- Table
  - gdcm::Table, [737](#)
- Table16
  - vtkLookupTable16, [911](#)
- TableEntry
  - gdcm::TableEntry, [738](#)
- TableReader
  - gdcm::TableReader, [739](#)
- TableRow
  - gdcm::network::TableRow, [741](#)
- Tag
  - gdcm::Tag, [743](#)
- tag
  - gdcm::Tag, [748](#)
- TagField
  - gdcm::DataElement, [291](#)
- TagMismatchError
  - gdcm::Parser, [555](#)
- TagPath
  - gdcm::TagPath, [749](#)
- TagToValue
  - gdcm::Scanner, [645](#)
- TagToValueValueType
  - gdcm::Scanner, [645](#)
- tags
  - gdcm::Tag, [748](#)
- TalairachBrainAtlasFrameofReference
  - gdcm::UIDs, [773](#)
- TestAbortOff
  - gdcm::SimpleSubjectWatcher, [685](#)
- TestAbortOn
  - gdcm::SimpleSubjectWatcher, [685](#)
- TestPBKDF2
  - gdcm::ASN1, [168](#)
- Testing
  - gdcm::Testing, [751](#)
- TestsList.txt, [1192](#)
- TextSRStorageTrialRetired
  - gdcm::UIDs, [776](#)
- ThreadedExecute
  - vtkImageRGBToYBR, [907](#)
  - vtkImageYBRToRGB, [909](#)
- ThreadedRequestData
  - vtkGDCMThreadedImageReader2, [891](#)
  - vtkImageMapToColors16, [901](#)
  - vtkImageMapToWindowLevelColors2, [904](#)
- to\_string
  - gdcm, [130](#)
- ToPyObject
  - gdcm::PythonFilter, [613](#)
- ToString
  - gdcm::StringFilter, [712](#)
- ToStringPair
  - gdcm::StringFilter, [712](#)
- ToUnixSlashes
  - gdcm::Filename, [386](#)
- ToWindowsSlashes
  - gdcm::Filename, [386](#)
- ToshibaPrivateDataStorage
  - gdcm::MediaStorage, [509](#)
- Trace
  - gdcm::Trace, [755](#)
- TransferSyntax
  - gdcm::TransferSyntax, [760](#)
- TransferSyntaxArrayType
  - gdcm::PresentationContext, [593](#)
- TransferSyntaxStringsType
  - gdcm::UIDs, [772](#)
- TransferSyntaxSub

- gdcm::network::TransferSyntaxSub, 762
- Transition
  - gdcm::network::Transition, 763
- transitions
  - gdcm::network::TableRow, 741
- Trim
  - gdcm::String, 710
- TrimInternal
  - gdcm::CodeString, 248
- Truncate
  - gdcm::String, 710
- TryJPEG2000Codec
  - gdcm::Bitmap, 216
  - gdcm::ImageChangeTransferSyntax, 428
- TryJPEG2000Codec2
  - gdcm::Bitmap, 216
- TryJPEGCodec
  - gdcm::Bitmap, 216
  - gdcm::ImageChangeTransferSyntax, 428
- TryJPEGCodec2
  - gdcm::Bitmap, 216
- TryJPEGLSCodec
  - gdcm::Bitmap, 216
  - gdcm::ImageChangeTransferSyntax, 428
- TryKAKADUCodec
  - gdcm::Bitmap, 216
- TryPVRGCodec
  - gdcm::Bitmap, 216
- TryRAWCodec
  - gdcm::Bitmap, 216
  - gdcm::ImageChangeTransferSyntax, 428
- TryRLECodec
  - gdcm::Bitmap, 216
  - gdcm::ImageChangeTransferSyntax, 428
- Type
  - gdcm::Element, 337
  - gdcm::Element< TVR, VM::VM1\_n >, 340
  - gdcm::Type, 765
  - gdcm::VL, 845
- TypeType
  - gdcm::Type, 765
- UI
  - gdcm::VR, 853
- UIComp
  - gdcm, 124
- UIDGenerator
  - gdcm::UIDGenerator, 767
- UINT12
  - gdcm::PixelFormat, 573
- UINT16
  - gdcm::PixelFormat, 574
- UINT32
  - gdcm::PixelFormat, 574
- UINT64
  - gdcm::PixelFormat, 574
- UINT8
  - gdcm::PixelFormat, 573
- UL
  - gdcm::VR, 853
- ULAction
  - gdcm::network::ULAction, 788
- ULActionAE6
  - gdcm::network::ULConnection, 820
- ULBasicCallback
  - gdcm::network::ULBasicCallback, 818
- ULConnection
  - gdcm::network::ULConnection, 819
- ULConnectionCallback
  - gdcm::network::ULConnectionCallback, 822
- ULConnectionInfo
  - gdcm::network::ULConnectionInfo, 823
- ULConnectionManager
  - gdcm::network::ULConnection, 821
  - gdcm::network::ULConnectionManager, 825
- ULError
  - gdcm::network::ULError, 827
- ULTransitionTable
  - gdcm::network::ULTransitionTable, 828
- ULWritingCallback
  - gdcm::network::ULWritingCallback, 829
- UN
  - gdcm::VR, 853
- UNKNOWN
  - gdcm::PhotometricInterpretation, 570
- UNKNOWN
  - gdcm::CSAHeader, 270
  - gdcm::LookupTable, 498
  - gdcm::Orientation, 545
  - gdcm::PixelFormat, 574
  - gdcm::Spacing, 695
  - gdcm::Surface, 719
  - gdcm::Type, 765
- URI
  - gdcm::MediaStorage, 510
- US
  - gdcm::VR, 853
- US\_SS
  - gdcm::VR, 853
- US\_SS\_OW
  - gdcm::VR, 853
- UT
  - gdcm::VR, 853
- UTComp
  - gdcm, 124
- uid\_1\_2\_840\_10008\_15\_0\_3\_1
  - gdcm::UIDs, 784
- uid\_1\_2\_840\_10008\_15\_0\_3\_10

gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_11  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_12  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_13  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_14  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_15  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_16  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_17  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_18  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_19  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_20  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_21  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_22  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_23  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_24  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_25  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_26  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_27  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_28  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_29  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_30  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_31  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_4  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_5  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_6  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_7

gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_8  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_3\_9  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_1  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_2  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_3  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_4  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_5  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_6  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_7  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_15\_0\_4\_8  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_1\_1  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_20\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_1\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_2  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_20\_2\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_1  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_1\_99  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_2  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_100  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_50  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_51  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_52  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_53  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_54  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_55  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_56



gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_57  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_58  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_59  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_60  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_61  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_62  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_63  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_64  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_65  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_66  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_70  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_80  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_81  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_90  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_91  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_92  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_93  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_94  
gdcm::UIDs, [779](#)  
uid\_1\_2\_840\_10008\_1\_2\_4\_95  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_5  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_6\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_2\_6\_2  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_3\_10  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_40  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_40\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_42  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_42\_1

gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_10  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_11  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_12  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_13  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_14  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_15  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_16  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_17  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_18  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_2  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_3  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_4  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_5  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_6  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_7  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_8  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_1\_9  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_2\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_4\_2\_2  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_1\_9  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_2\_16\_4  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_2\_6\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_1\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_1\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_1\_4  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_2\_1



gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_3\_1  
gdcm::UIDs, [780](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_3\_2  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_3\_3  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_3\_4  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_3\_5  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_5\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_5\_4  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_5\_5  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_3\_1\_2\_6\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_4\_2  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_14  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_15  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_16  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_16\_376  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_17  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_17\_376  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_18  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_18\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_2  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_22  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_23  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_24  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_24\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_25  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_26  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_27

gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_29  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_30  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_31  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_32  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_33  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_4  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_4\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_4\_2  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_9  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_1\_9\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_10  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_104\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_3  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_11\_4  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_128  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_129  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_1\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_2\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_12\_3  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_1

gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_13\_1\_3  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_1\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_2\_1  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3  
gdcm::UIDs, [781](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_1\_3\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_20  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_2\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_3\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_3  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_4  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_5  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_6  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_7  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_8  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_481\_9  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_4\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_5  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6

gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_3  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_4  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_66\_5  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_67  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_6\_2  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_2\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_3  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_4\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_3  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_5\_4  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_1\_6  
gdcm::UIDs, [785](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_77\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_3

gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_7\_4  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_8  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_11  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_2  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_22  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_3  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_33  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_4  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_40  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_50  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_59  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_65  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_88\_67  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_2  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_1\_3  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_2\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_3\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_1\_9\_4\_1  
gdcm::UIDs, [782](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_1  
gdcm::UIDs, [783](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_1\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_2

gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_2\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_1\_2\_3\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_31  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_32  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_32\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_32\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_32\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_33  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_4\_4  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_34\_5  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_37\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_37\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_37\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_38\_1  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_38\_2  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_38\_3  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_41  
gdcm::UIDs, [784](#)  
uid\_1\_2\_840\_10008\_5\_1\_4\_42

- gdcmm::UIDs, [784](#)
- UltrasoundImageStorage
  - gdcmm::MediaStorage, [508](#)
  - gdcmm::UIDs, [775](#)
- UltrasoundImageStorageRetired
  - gdcmm::MediaStorage, [508](#)
  - gdcmm::UIDs, [775](#)
- UltrasoundMultiFrameImageStorage
  - gdcmm::MediaStorage, [508](#)
- UltrasoundMultiFrameImageStorageRetired
  - gdcmm::MediaStorage, [508](#)
- UltrasoundMultiframeImageStorage
  - gdcmm::UIDs, [775](#)
- UltrasoundMultiframeImageStorageRetired
  - gdcmm::UIDs, [775](#)
- UnInstallPipeline
  - vtkImageColorViewer, [898](#)
- UnRegister
  - gdcmm::Object, [536](#)
- UndefinedEntityError
  - gdcmm::Parser, [555](#)
- underline
  - gdcmm::terminal, [137](#)
- UnexpectedStateError
  - gdcmm::Parser, [555](#)
- UnifiedProcedureStepEventSOPClass
  - gdcmm::UIDs, [777](#)
- UnifiedProcedureStepPullSOPClass
  - gdcmm::UIDs, [777](#)
- UnifiedProcedureStepPushSOPClass
  - gdcmm::UIDs, [777](#)
- UnifiedProcedureStepWatchSOPClass
  - gdcmm::UIDs, [777](#)
- UnifiedWorklistandProcedureStepSOPInstance
  - gdcmm::UIDs, [777](#)
- UnifiedWorklistandProcedureStepServiceClass
  - gdcmm::UIDs, [777](#)
- Unknown
  - gdcmm::SwapCode, [730](#)
  - gdcmm::TransferSyntax, [759](#)
- Unpack
  - gdcmm::Unpacker12Bits, [834](#)
- Update
  - gdcmm::Curve, [282](#)
  - gdcmm::Overlay, [552](#)
- UpdateDisplayExtent
  - vtkImageColorViewer, [898](#)
- UpdateOrientation
  - vtkImageColorViewer, [898](#)
- UpdatePosition
  - gdcmm::ByteBuffer, [222](#)
- Usage
  - gdcmm::Usage, [835](#)
- UsageType
  - gdcmm::Usage, [835](#)
- UseDictAlways
  - gdcmm::PythonFilter, [613](#)
  - gdcmm::StringFilter, [713](#)
- UserInformation
  - gdcmm::network::UserInformation, [837](#)
- UserOption
  - gdcmm::Usage, [835](#)
- UserOrdering
  - gdcmm::SerieHelper, [672](#)
- V
  - gdcmm::Validate, [840](#)
- VERBOSE\_STYLE
  - gdcmm::Printer, [603](#)
- VERTEX
  - gdcmm::MeshPrimitive, [518](#)
- VIEWType
  - gdcmm::Surface, [719](#)
- VIEWType\_END
  - gdcmm::Surface, [719](#)
- VL
  - gdcmm::VL, [845](#)
- VL16
  - gdcmm::VR, [853](#)
- VL32
  - gdcmm::VR, [853](#)
- VLEndoscopicImageStorage
  - gdcmm::MediaStorage, [510](#)
  - gdcmm::UIDs, [776](#)
- VLImageStorageTrialRetired
  - gdcmm::UIDs, [776](#)
- VLMicroscopicImageStorage
  - gdcmm::MediaStorage, [510](#)
  - gdcmm::UIDs, [776](#)
- VLMultiframeImageStorageTrialRetired
  - gdcmm::UIDs, [776](#)
- VLPhotographicImageStorage
  - gdcmm::MediaStorage, [510](#)
  - gdcmm::UIDs, [776](#)
- VLSlideCoordinatesMicroscopicImageStorage
  - gdcmm::UIDs, [776](#)
- VLWholeSlideMicroscopyImageStorage
  - gdcmm::MediaStorage, [510](#)
  - gdcmm::UIDs, [779](#)
- VM
  - gdcmm::VM, [849](#)
- VM0
  - gdcmm::VM, [848](#)
- VM1
  - gdcmm::VM, [848](#)
- VM10
  - gdcmm::VM, [848](#)
- VM12

- gdcM::VM, [848](#)
- VM16
  - gdcM::VM, [848](#)
- VM18
  - gdcM::VM, [848](#)
- VM1\_2
  - gdcM::VM, [849](#)
- VM1\_3
  - gdcM::VM, [849](#)
- VM1\_32
  - gdcM::VM, [849](#)
- VM1\_4
  - gdcM::VM, [849](#)
- VM1\_5
  - gdcM::VM, [849](#)
- VM1\_8
  - gdcM::VM, [849](#)
- VM1\_99
  - gdcM::VM, [849](#)
- VM1\_n
  - gdcM::VM, [849](#)
- VM2
  - gdcM::VM, [848](#)
- VM24
  - gdcM::VM, [848](#)
- VM256
  - gdcM::VM, [849](#)
- VM28
  - gdcM::VM, [848](#)
- VM2\_2n
  - gdcM::VM, [849](#)
- VM2\_n
  - gdcM::VM, [849](#)
- VM3
  - gdcM::VM, [848](#)
- VM30\_30n
  - gdcM::VM, [849](#)
- VM32
  - gdcM::VM, [848](#)
- VM35
  - gdcM::VM, [848](#)
- VM3\_3n
  - gdcM::VM, [849](#)
- VM3\_4
  - gdcM::VM, [849](#)
- VM3\_n
  - gdcM::VM, [849](#)
- VM4
  - gdcM::VM, [848](#)
- VM47\_47n
  - gdcM::VM, [849](#)
- VM4\_4n
  - gdcM::VM, [849](#)
- VM5
  - gdcM::VM, [848](#)
- VM6
  - gdcM::VM, [848](#)
- VM6\_6n
  - gdcM::VM, [849](#)
- VM7\_7n
  - gdcM::VM, [849](#)
- VM8
  - gdcM::VM, [848](#)
- VM9
  - gdcM::VM, [848](#)
- VM99
  - gdcM::VM, [849](#)
- VM\_END
  - gdcM::VM, [849](#)
- VMType
  - gdcM::Attribute, [171](#)
  - gdcM::Attribute< Group, Element, TVR, VM::VM1 >, [177](#)
  - gdcM::VM, [848](#)
- VOILUTBoxSOPClass
  - gdcM::UIDs, [775](#)
- VR
  - gdcM::VR, [854](#)
- VR\_END
  - gdcM::VR, [853](#)
- VR\_VM1
  - gdcM::VR, [853](#)
- VRALL
  - gdcM::VR, [853](#)
- VRASCII
  - gdcM::VR, [853](#)
- VRBINARY
  - gdcM, [131](#)
  - gdcM::VR, [853](#)
- VRField
  - gdcM::CSAElement, [268](#)
  - gdcM::DataElement, [291](#)
- VRType
  - gdcM::VR, [852](#)
- VRTypeTemplateCase
  - gdcMVR.h, [1187](#)
- VT100
  - gdcM::terminal, [138](#)
- VTK\_CMYK
  - vtkGDCMImageReader.h, [1193](#)
  - vtkGDCMImageReader2.h, [1194](#)
- VTK\_INVERSE\_LUMINANCE
  - vtkGDCMImageReader.h, [1193](#)
  - vtkGDCMImageReader2.h, [1194](#)
- VTK\_LEGACY
  - vtkImageColorViewer, [898](#)
- VTK\_LOOKUP\_TABLE
  - vtkGDCMImageReader.h, [1193](#)

- vtkGDCMImageReader2.h, [1194](#)
- VTK\_YBR
  - vtkGDCMImageReader.h, [1193](#)
  - vtkGDCMImageReader2.h, [1194](#)
- Valid
  - gdcm::Preamble, [592](#)
- Validate
  - gdcm::PixelFormat, [576](#)
  - gdcm::Validate, [839](#)
- ValidateQuery
  - gdcm::BaseRootQuery, [203](#)
  - gdcm::FindPatientRootQuery, [399](#)
  - gdcm::FindStudyRootQuery, [401](#)
  - gdcm::MovePatientRootQuery, [528](#)
  - gdcm::MoveStudyRootQuery, [530](#)
- Validation
  - gdcm::Validate, [840](#)
- Value
  - gdcm::Value, [841](#)
- value
  - gdcm::STATIC\_ASSERTION\_FAILURE< true >, [699](#)
  - gdcm::SerieHelper::Rule, [642](#)
- value\_type
  - gdcm::CodeString, [247](#)
  - gdcm::LO, [495](#)
  - gdcm::String, [709](#)
- ValueField
  - gdcm::DataElement, [291](#)
  - gdcm::PDBelement, [560](#)
- ValueLengthField
  - gdcm::DataElement, [291](#)
- ValueMultiplicityField
  - gdcm::CSAElement, [268](#)
- ValuePtr
  - gdcm::DataElement, [285](#)
- ValueType
  - gdcm::Scanner, [645](#)
- VerificationSOPClass
  - gdcm::UIDs, [772](#)
- Verify
  - gdcm::Defs, [308](#)
  - gdcm::Macro, [502](#)
  - gdcm::Module, [522](#)
- Version
  - gdcm::Version, [843](#)
- Video
  - gdcm::MediaStorage, [510](#)
- VideoEndoscopicImageStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [776](#)
- VideoMicroscopicImageStorage
  - gdcm::UIDs, [776](#)
- VideoPhotographicImageStorage
  - gdcm::UIDs, [776](#)
- vtkBooleanMacro
  - vtkGDCMImageReader, [863](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMThreadedImageReader, [888](#)
  - vtkGDCMThreadedImageReader2, [891](#)
  - vtkImageColorViewer, [898](#)
  - vtkImageMapToColors16, [901](#)
- vtkGDCMImageReader, [859](#)
  - ~vtkGDCMImageReader, [861](#)
  - ApplyInverseVideo, [864](#)
  - ApplyLookupTable, [864](#)
  - ApplyPlanarConfiguration, [864](#)
  - ApplyShiftScale, [864](#)
  - ApplyYBRToRGB, [864](#)
  - CanReadFile, [862](#)
  - Curve, [864](#)
  - DirectionCosines, [864](#)
  - ExecuteData, [862](#)
  - ExecuteInformation, [862](#)
  - FileNames, [864](#)
  - FillMedicalImageInformation, [862](#)
  - ForceRescale, [864](#)
  - GetDescriptiveName, [862](#)
  - GetFileExtensions, [862](#)
  - GetIconImage, [862](#)
  - GetOverlay, [862](#)
  - IconDataScalarType, [864](#)
  - IconImageDataExtent, [864](#)
  - IconNumberOfScalarComponents, [864](#)
  - ImageFormat, [864](#)
  - ImageOrientationPatient, [864](#)
  - ImagePositionPatient, [864](#)
  - LoadIconImage, [865](#)
  - LoadOverlays, [865](#)
  - LoadSingleFile, [862](#)
  - LossyFlag, [865](#)
  - MedicalImageProperties, [865](#)
  - New, [862](#)
  - NumberOfIconImages, [865](#)
  - NumberOfOverlays, [865](#)
  - PlanarConfiguration, [865](#)
  - PrintSelf, [862](#)
  - RequestDataCompat, [862](#)
  - RequestInformationCompat, [862](#)
  - Scale, [865](#)
  - SetCurve, [862](#)
  - SetFileNames, [862](#)
  - SetFilePattern, [863](#)
  - SetFilePrefix, [863](#)
  - SetMedicalImageProperties, [863](#)
  - Shift, [865](#)
  - vtkBooleanMacro, [863](#)

- vtkGDCMImageReader, 861
- vtkGDCMMedicalImageProperties, 877
- vtkGetMacro, 863
- vtkGetObjectMacro, 863
- vtkGetStringMacro, 864
- vtkGetVector3Macro, 864
- vtkGetVector6Macro, 864
- vtkSetMacro, 864
- vtkSetVector6Macro, 864
- vtkTypeRevisionMacro, 864
- vtkGDCMImageReader.h, 1192
  - VTK\_CMYK, 1193
  - VTK\_INVERSE\_LUMINANCE, 1193
  - VTK\_LOOKUP\_TABLE, 1193
  - VTK\_YBR, 1193
- vtkGDCMImageReader2, 865
  - ~vtkGDCMImageReader2, 868
  - ApplyInverseVideo, 870
  - ApplyLookupTable, 870
  - ApplyPlanarConfiguration, 870
  - ApplyShiftScale, 870
  - ApplyYBRToRGB, 870
  - CanReadFile, 868
  - Curve, 870
  - DirectionCosines, 870
  - FillMedicalImageInformation, 868
  - ForceRescale, 870
  - GetDescriptiveName, 868
  - GetFileExtensions, 868
  - GetIconImage, 868
  - GetIconImagePort, 868
  - GetOverlay, 868
  - GetOverlayPort, 868
  - IconDataScalarType, 870
  - IconImageDataExtent, 870
  - IconNumberOfScalarComponents, 870
  - ImageFormat, 870
  - ImageOrientationPatient, 870
  - ImagePositionPatient, 870
  - LoadIconImage, 870
  - LoadOverlays, 870
  - LoadSingleFile, 868
  - LossyFlag, 870
  - New, 868
  - NumberOfIconImages, 870
  - NumberOfOverlays, 870
  - PlanarConfiguration, 871
  - PrintSelf, 868
  - ProcessRequest, 868
  - RequestData, 868
  - RequestDataCompat, 868
  - RequestInformation, 868
  - RequestInformationCompat, 868
  - Scale, 871
  - SetCurve, 869
  - SetFilePattern, 869
  - SetFilePrefix, 869
  - SetMedicalImageProperties, 869
  - Shift, 871
  - vtkBooleanMacro, 869
  - vtkGDCMImageReader2, 868
  - vtkGDCMMedicalImageProperties, 877
  - vtkGetMacro, 869
  - vtkGetObjectMacro, 869
  - vtkGetStringMacro, 869
  - vtkGetVector3Macro, 869
  - vtkGetVector6Macro, 869
  - vtkSetMacro, 870
  - vtkSetVector6Macro, 870
  - vtkTypeRevisionMacro, 870
- vtkGDCMImageReader2.h, 1193
  - VTK\_CMYK, 1194
  - VTK\_INVERSE\_LUMINANCE, 1194
  - VTK\_LOOKUP\_TABLE, 1194
  - VTK\_YBR, 1194
- vtkGDCMImageWriter, 871
  - ~vtkGDCMImageWriter, 873
  - CompressionTypes, 873
  - GetDescriptiveName, 873
  - GetFileExtensions, 873
  - GetFileName, 873
  - JPEG2000\_COMPRESSION, 873
  - JPEG\_COMPRESSION, 873
  - JPEGLS\_COMPRESSION, 873
  - NO\_COMPRESSION, 873
  - New, 873
  - PrintSelf, 873
  - RLE\_COMPRESSION, 873
  - SetDirectionCosines, 874
  - SetDirectionCosinesFromImageOrientationPatient, 874
  - SetFileNames, 874
  - SetMedicalImageProperties, 874
  - vtkBooleanMacro, 874
  - vtkGDCMImageWriter, 873
  - vtkGDCMMedicalImageProperties, 877
  - vtkGetMacro, 874
  - vtkGetObjectMacro, 874
  - vtkGetStringMacro, 874
  - vtkSetMacro, 875
  - vtkSetStringMacro, 875
  - vtkTypeRevisionMacro, 875
  - Write, 875
  - WriteGDCMData, 875
  - WriteSlice, 875
- vtkGDCMImageWriter.h, 1194
- vtkGDCMMedicalImageProperties, 875
  - ~vtkGDCMMedicalImageProperties, 877



- Clear, [877](#)
- GetFile, [877](#)
- New, [877](#)
- PrintSelf, [877](#)
- PushBackFile, [877](#)
- vtkGDCMImageReader, [877](#)
- vtkGDCMImageReader2, [877](#)
- vtkGDCMImageWriter, [877](#)
- vtkGDCMMedicalImageProperties, [877](#)
- vtkTypeRevisionMacro, [877](#)
- vtkGDCMMedicalImageProperties.h, [1194](#)
- vtkGDCMPolyDataReader, [877](#)
  - ~vtkGDCMPolyDataReader, [879](#)
  - FileName, [880](#)
  - FillMedicalImageInformation, [879](#)
  - MedicalImageProperties, [880](#)
  - New, [879](#)
  - PrintSelf, [879](#)
  - RTStructSetProperties, [880](#)
  - RequestData, [879](#)
  - RequestData\_HemodynamicWaveformStorage, [879](#)
  - RequestData\_RTStructureSetStorage, [879](#)
  - RequestInformation, [880](#)
  - RequestInformation\_HemodynamicWaveformStorage, [880](#)
  - RequestInformation\_RTStructureSetStorage, [880](#)
  - vtkGDCMPolyDataReader, [879](#)
  - vtkGetObjectMacro, [880](#)
  - vtkGetStringMacro, [880](#)
  - vtkSetStringMacro, [880](#)
  - vtkTypeRevisionMacro, [880](#)
- vtkGDCMPolyDataReader.h, [1195](#)
- vtkGDCMPolyDataWriter, [880](#)
  - ~vtkGDCMPolyDataWriter, [882](#)
  - InitializeRTStructSet, [882](#)
  - MedicalImageProperties, [883](#)
  - New, [882](#)
  - PrintSelf, [882](#)
  - RTStructSetProperties, [883](#)
  - SetMedicalImageProperties, [882](#)
  - SetNumberOfInputPorts, [882](#)
  - SetRTStructSetProperties, [883](#)
  - vtkGDCMPolyDataWriter, [882](#)
  - vtkTypeRevisionMacro, [883](#)
  - WriteData, [883](#)
  - WriteRTSTRUCTData, [883](#)
  - WriteRTSTRUCTInfo, [883](#)
- vtkGDCMPolyDataWriter.h, [1196](#)
- vtkGDCMTesting, [883](#)
  - ~vtkGDCMTesting, [885](#)
  - GetGDCMDataRoot, [885](#)
  - GetMD5MetalImage, [885](#)
  - GetMHDMD5FromFile, [885](#)
  - GetNumberOfMD5MetalImages, [885](#)
  - GetRAWMD5FromFile, [885](#)
  - GetVTKDataRoot, [885](#)
  - MD5MetalImagesType, [885](#)
  - New, [885](#)
  - PrintSelf, [885](#)
  - vtkGDCMTesting, [885](#)
  - vtkTypeRevisionMacro, [886](#)
- vtkGDCMTesting.h, [1196](#)
- vtkGDCMThreadedImageReader, [886](#)
  - ~vtkGDCMThreadedImageReader, [887](#)
  - ExecuteData, [888](#)
  - ExecuteInformation, [888](#)
  - New, [888](#)
  - PrintSelf, [888](#)
  - ReadFiles, [888](#)
  - RequestDataCompat, [888](#)
  - vtkBooleanMacro, [888](#)
  - vtkGDCMThreadedImageReader, [887](#)
  - vtkGetMacro, [888](#)
  - vtkSetMacro, [888](#)
  - vtkTypeRevisionMacro, [888](#)
- vtkGDCMThreadedImageReader.h, [1197](#)
- vtkGDCMThreadedImageReader2, [888](#)
  - ~vtkGDCMThreadedImageReader2, [890](#)
  - GetFileName, [890](#)
  - New, [890](#)
  - PrintSelf, [890](#)
  - RequestInformation, [890](#)
  - SetFileName, [890](#)
  - SetFileNames, [890](#)
  - SplitExtent, [890](#)
  - ThreadedRequestData, [891](#)
  - vtkBooleanMacro, [891](#)
  - vtkGDCMThreadedImageReader2, [890](#)
  - vtkGetMacro, [891](#)
  - vtkGetObjectMacro, [891](#)
  - vtkGetVector3Macro, [891](#)
  - vtkGetVector6Macro, [891](#)
  - vtkSetMacro, [891](#)
  - vtkSetVector3Macro, [891](#)
  - vtkSetVector6Macro, [892](#)
  - vtkTypeRevisionMacro, [892](#)
- vtkGDCMThreadedImageReader2.h, [1198](#)
- vtkGetMacro
  - vtkGDCMImageReader, [863](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMThreadedImageReader, [888](#)
  - vtkGDCMThreadedImageReader2, [891](#)
  - vtkImageColorViewer, [898](#)
  - vtkImageMapToColors16, [901](#)
  - vtkImageMapToWindowLevelColors2, [904](#)
- vtkGetObjectMacro
  - vtkGDCMImageReader, [863](#)



- vtkGDCMImageReader2, [869](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMPolyDataReader, [880](#)
  - vtkGDCMThreadedImageReader2, [891](#)
  - vtkImageColorViewer, [898](#)
  - vtkImageMapToColors16, [901](#)
- vtkGetStringMacro
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMImageWriter, [874](#)
  - vtkGDCMPolyDataReader, [880](#)
  - vtkRTStructSetProperties, [915](#)
- vtkGetVector3Macro
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMThreadedImageReader2, [891](#)
- vtkGetVector6Macro
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [869](#)
  - vtkGDCMThreadedImageReader2, [891](#)
- vtkImageColorViewer, [892](#)
  - ~vtkImageColorViewer, [895](#)
  - AddInput, [895](#)
  - AddInputConnection, [895](#)
  - FirstRender, [898](#)
  - GetColorLevel, [895](#)
  - GetColorWindow, [895](#)
  - GetInput, [895](#)
  - GetOffScreenRendering, [895](#)
  - GetOverlayVisibility, [896](#)
  - GetPosition, [896](#)
  - GetSize, [896](#)
  - GetSliceMax, [896](#)
  - GetSliceMin, [896](#)
  - GetSliceRange, [896](#)
  - GetWindowName, [896](#)
  - ImageActor, [898](#)
  - InstallPipeline, [896](#)
  - Interactor, [898](#)
  - InteractorStyle, [898](#)
  - New, [896](#)
  - OverlayImageActor, [899](#)
  - PrintSelf, [896](#)
  - Render, [896](#)
  - RenderWindow, [899](#)
  - Renderer, [899](#)
  - SLICE\_ORIENTATION\_XY, [895](#)
  - SLICE\_ORIENTATION\_XZ, [895](#)
  - SLICE\_ORIENTATION\_YZ, [895](#)
  - SetColorLevel, [896](#)
  - SetColorWindow, [896](#)
  - SetDisplayId, [896](#)
  - SetInput, [896](#)
  - SetInputConnection, [896](#)
  - SetOffScreenRendering, [897](#)
  - SetOverlayVisibility, [897](#)
  - SetParentId, [897](#)
  - SetPosition, [897](#)
  - SetRenderWindow, [897](#)
  - SetRenderer, [897](#)
  - SetSize, [897](#)
  - SetSlice, [897](#)
  - SetSliceOrientation, [897](#)
  - SetSliceOrientationToXY, [897](#)
  - SetSliceOrientationToXZ, [897](#)
  - SetSliceOrientationToYZ, [897](#)
  - SetWindowId, [898](#)
  - SetupInteractor, [897](#)
  - Slice, [899](#)
  - SliceOrientation, [899](#)
  - UnInstallPipeline, [898](#)
  - UpdateDisplayExtent, [898](#)
  - UpdateOrientation, [898](#)
  - VTK\_LEGACY, [898](#)
  - vtkBooleanMacro, [898](#)
  - vtkGetMacro, [898](#)
  - vtkGetObjectMacro, [898](#)
  - vtkImageColorViewer, [895](#)
  - vtkImageColorViewerCallback, [898](#)
  - vtkTypeRevisionMacro, [898](#)
  - WindowLevel, [899](#)
- vtkImageColorViewer.h, [1198](#)
- vtkImageColorViewerCallback
  - vtkImageColorViewer, [898](#)
- vtkImageMapToColors16, [899](#)
  - ~vtkImageMapToColors16, [901](#)
  - ActiveComponent, [902](#)
  - DataWasPassed, [902](#)
  - GetMTime, [901](#)
  - LookupTable, [902](#)
  - New, [901](#)
  - OutputFormat, [902](#)
  - PassAlphaToOutput, [902](#)
  - PrintSelf, [901](#)
  - RequestData, [901](#)
  - RequestInformation, [901](#)
  - SetLookupTable, [901](#)
  - SetOutputFormatToLuminance, [901](#)
  - SetOutputFormatToLuminanceAlpha, [901](#)
  - SetOutputFormatToRGB, [901](#)
  - SetOutputFormatToRGBA, [901](#)
  - ThreadedRequestData, [901](#)
  - vtkBooleanMacro, [901](#)
  - vtkGetMacro, [901](#)
  - vtkGetObjectMacro, [901](#)
  - vtkImageMapToColors16, [901](#)
  - vtkSetMacro, [901](#), [902](#)
  - vtkTypeRevisionMacro, [902](#)

- vtkImageMapToColors16.h, 1199
- vtkImageMapToWindowLevelColors2, 902
  - ~vtkImageMapToWindowLevelColors2, 904
  - Level, 904
  - New, 904
  - PrintSelf, 904
  - RequestData, 904
  - RequestInformation, 904
  - ThreadedRequestData, 904
  - vtkGetMacro, 904
  - vtkImageMapToWindowLevelColors2, 904
  - vtkSetMacro, 904
  - vtkTypeRevisionMacro, 904
  - Window, 904
- vtkImageMapToWindowLevelColors2.h, 1199
- vtkImagePlanarComponentsToComponents, 904
  - ~vtkImagePlanarComponentsToComponents, 906
  - New, 906
  - PrintSelf, 906
  - RequestData, 906
  - vtkImagePlanarComponentsToComponents, 906
  - vtkTypeRevisionMacro, 906
- vtkImagePlanarComponentsToComponents.h, 1200
- vtkImageRGBToYBR, 906
  - ~vtkImageRGBToYBR, 907
  - New, 907
  - PrintSelf, 907
  - ThreadedExecute, 907
  - vtkImageRGBToYBR, 907
  - vtkTypeRevisionMacro, 907
- vtkImageRGBToYBR.h, 1200
- vtkImageYBRToRGB, 908
  - ~vtkImageYBRToRGB, 909
  - New, 909
  - PrintSelf, 909
  - ThreadedExecute, 909
  - vtkImageYBRToRGB, 909
  - vtkTypeRevisionMacro, 909
- vtkImageYBRToRGB.h, 1201
- vtkLookupTable16, 909
  - ~vtkLookupTable16, 910
  - Build, 910
  - GetPointer, 911
  - MapScalarsThroughTable2, 911
  - New, 911
  - PrintSelf, 911
  - SetNumberOfTableValues, 911
  - Table16, 911
  - vtkLookupTable16, 910
  - vtkTypeRevisionMacro, 911
  - WritePointer, 911
- vtkLookupTable16.h, 1201
- vtkRTStructSetProperties, 911
  - ~vtkRTStructSetProperties, 913
  - AddContourReferencedFrameOfReference, 913
  - AddReferencedFrameOfReference, 914
  - AddStructureSetROI, 914
  - AddStructureSetROIObservation, 914
  - Clear, 914
  - DeepCopy, 914
  - GetContourReferencedFrameOfReferenceClassUID, 914
  - GetContourReferencedFrameOfReferenceInstanceUID, 914
  - GetNumberOfContourReferencedFrameOfReferences, 914
  - GetNumberOfReferencedFrameOfReferences, 914
  - GetNumberOfStructureSetROIs, 914
  - GetReferencedFrameOfReferenceClassUID, 914
  - GetReferencedFrameOfReferenceInstanceUID, 914
  - GetStructureSetObservationNumber, 914
  - GetStructureSetROIDescription, 914
  - GetStructureSetROIGenerationAlgorithm, 914
  - GetStructureSetROIName, 914
  - GetStructureSetROINumber, 914
  - GetStructureSetROIObservationLabel, 914
  - GetStructureSetROIRefFrameRefUID, 914
  - GetStructureSetRTROIInterpretedType, 914
  - Internals, 915
  - New, 914
  - PrintSelf, 915
  - ReferenceFrameOfReferenceUID, 915
  - ReferenceSeriesInstanceUID, 915
  - SOPInstanceUID, 916
  - SeriesInstanceUID, 916
  - StructureSetDate, 916
  - StructureSetLabel, 916
  - StructureSetName, 916
  - StructureSetTime, 916
  - StudyInstanceUID, 916
  - vtkGetStringMacro, 915
  - vtkRTStructSetProperties, 913
  - vtkSetStringMacro, 915
  - vtkTypeRevisionMacro, 915
- vtkRTStructSetProperties.h, 1202
- vtkSetMacro
  - vtkGDCMImageReader, 864
  - vtkGDCMImageReader2, 870
  - vtkGDCMImageWriter, 875
  - vtkGDCMThreadedImageReader, 888
  - vtkGDCMThreadedImageReader2, 891
  - vtkImageMapToColors16, 901, 902
  - vtkImageMapToWindowLevelColors2, 904
- vtkSetStringMacro
  - vtkGDCMImageWriter, 875
  - vtkGDCMPolyDataReader, 880
  - vtkRTStructSetProperties, 915
- vtkSetVector3Macro

- vtkGDCMThreadedImageReader2, [891](#)
- vtkSetVector6Macro
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
  - vtkGDCMThreadedImageReader2, [892](#)
- vtkTypeRevisionMacro
  - vtkGDCMImageReader, [864](#)
  - vtkGDCMImageReader2, [870](#)
  - vtkGDCMImageWriter, [875](#)
  - vtkGDCMMedicalImageProperties, [877](#)
  - vtkGDCMPolyDataReader, [880](#)
  - vtkGDCMPolyDataWriter, [883](#)
  - vtkGDCMTesting, [886](#)
  - vtkGDCMThreadedImageReader, [888](#)
  - vtkGDCMThreadedImageReader2, [892](#)
  - vtkImageColorViewer, [898](#)
  - vtkImageMapToColors16, [902](#)
  - vtkImageMapToWindowLevelColors2, [904](#)
  - vtkImagePlanarComponentsToComponents, [906](#)
  - vtkImageRGBToYBR, [907](#)
  - vtkImageYBRToRGB, [909](#)
  - vtkLookupTable16, [911](#)
  - vtkRTStructSetProperties, [915](#)
- WIREFRAME
  - gdcm::Surface, [719](#)
- WarningOff
  - gdcm::Trace, [757](#)
- WarningOn
  - gdcm::Trace, [757](#)
- Waveform
  - gdcm::MediaStorage, [510](#)
  - gdcm::Waveform, [916](#)
- WaveformStorageTrialRetired
  - gdcm::UIDs, [775](#)
- what
  - gdcm::Exception, [360](#)
- white
  - gdcm::terminal, [138](#)
- Window
  - vtkImageMapToWindowLevelColors2, [904](#)
- WindowLevel
  - vtkImageColorViewer, [899](#)
- Write
  - gdcm::ByteValue, [229](#)
  - gdcm::CSAHeader, [272](#)
  - gdcm::CommandDataSet, [252](#)
  - gdcm::DataElement, [290](#)
  - gdcm::DataSet, [301](#)
  - gdcm::Element, [337](#)
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
  - gdcm::EncodingImplementation< VR::VRASCII >, [354](#)
  - gdcm::EncodingImplementation< VR::VRBINARY >, [355](#)
  - gdcm::ExplicitDataElement, [363](#)
  - gdcm::File, [368](#)
  - gdcm::FileAnonymizer, [372](#)
  - gdcm::FileMetaInformation, [383](#)
  - gdcm::Fragment, [403](#)
  - gdcm::ImageWriter, [453](#)
  - gdcm::ImplicitDataElement, [457](#)
  - gdcm::Item, [470](#)
  - gdcm::PGXCodec, [569](#)
  - gdcm::PNMCodec, [590](#)
  - gdcm::PixmapWriter, [588](#)
  - gdcm::Preamble, [592](#)
  - gdcm::SegmentWriter, [658](#)
  - gdcm::SequenceOfFragments, [663](#)
  - gdcm::SequenceOfItems, [669](#)
  - gdcm::StreamImageWriter, [705](#)
  - gdcm::SurfaceWriter, [729](#)
  - gdcm::Tag, [747](#)
  - gdcm::VL, [846](#)
  - gdcm::VR, [855](#)
  - gdcm::VRVLSize< 0 >, [858](#)
  - gdcm::VRVLSize< 1 >, [858](#)
  - gdcm::ValueIO, [843](#)
  - gdcm::Writer, [921](#)
  - gdcm::network::AAAbortPDU, [141](#)
  - gdcm::network::AAAssociateACPDU, [143](#)
  - gdcm::network::AAAssociateRJPDU, [145](#)
  - gdcm::network::AAAssociateRQPDU, [149](#)
  - gdcm::network::AReleaseRPPDU, [164](#)
  - gdcm::network::AReleaseRQPDU, [166](#)
  - gdcm::network::AbstractSyntax, [151](#)
  - gdcm::network::ApplicationContext, [161](#)
  - gdcm::network::AsynchronousOperationsWindow< Sub, [169](#)
  - gdcm::network::BasePDU, [200](#)
  - gdcm::network::ImplementationClassUIDSub, [454](#)
  - gdcm::network::ImplementationUIDSub, [455](#)
  - gdcm::network::ImplementationVersionNameSub, [455](#)
  - gdcm::network::MaximumLengthSub, [504](#)
  - gdcm::network::PDataTFPDU, [558](#)
  - gdcm::network::PresentationContextAC, [595](#)
  - gdcm::network::PresentationContextRQ, [599](#)
  - gdcm::network::PresentationDataValue, [601](#)
  - gdcm::network::RoleSelectionSub, [641](#)
  - gdcm::network::SOPClassExtendedNegociationSub, [689](#)
  - gdcm::network::ServiceClassApplicationInformation, [674](#)
  - gdcm::network::TransferSyntaxSub, [762](#)
  - gdcm::network::UserInformation, [838](#)
  - vtkGDCMImageWriter, [875](#)

- Write16
  - gdcm::VL, [846](#)
- WriteASCII
  - gdcm::Element< TVR, VM::VM1\_n >, [341](#)
- WriteBuffer
  - gdcm::ByteValue, [229](#)
  - gdcm::SequenceOfFragments, [663](#)
- WriteBufferAsRGBA
  - gdcm::LookupTable, [500](#)
- WriteData
  - vtkGDCMPolyDataWriter, [883](#)
- WriteFooter
  - gdcm::DictConverter, [317](#)
- WriteGDCMData
  - vtkGDCMImageWriter, [875](#)
- WriteHeader
  - gdcm::DictConverter, [318](#)
- WriteHelpFile
  - gdcm::BaseRootQuery, [204](#)
- WriteImageInformation
  - gdcm::StreamImageWriter, [706](#)
- WriteImageSubregionRAW
  - gdcm::StreamImageWriter, [706](#)
- WritePointer
  - vtkLookupTable16, [911](#)
- WriteQuery
  - gdcm::BaseRootQuery, [204](#)
- WriteRTSTRUCTData
  - vtkGDCMPolyDataWriter, [883](#)
- WriteRTSTRUCTInfo
  - vtkGDCMPolyDataWriter, [883](#)
- WriteRawHeader
  - gdcm::StreamImageWriter, [706](#)
- WriteSlice
  - vtkGDCMImageWriter, [875](#)
- Writer
  - gdcm::Writer, [920](#)
- XML
  - gdcm::Printer, [603](#)
- XMLDictReader
  - gdcm::XMLDictReader, [923](#)
- XMLEncoding
  - gdcm::UIDs, [773](#)
- XMLPrinter
  - gdcm::XMLPrinter, [925](#)
- XMLPrivateDictReader
  - gdcm::XMLPrivateDictReader, [927](#)
- XRay3DAngiographicImageStorage
  - gdcm::MediaStorage, [510](#)
  - gdcm::UIDs, [776](#)
- XRay3DCraniofacialImageStorage
  - gdcm::UIDs, [776](#)
- XRayAngiographicBiPlaneImageStorageRetired
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [776](#)
- XRayAngiographicImageStorage
  - gdcm::MediaStorage, [509](#)
  - gdcm::UIDs, [776](#)
- XRayRadiationDoseSR
  - gdcm::MediaStorage, [510](#)
- XRayRadiationDoseSRStorage
  - gdcm::UIDs, [777](#)
- XRayRadiofluoroscopicImageStorage
  - gdcm::UIDs, [776](#)
- XRayRadiofluoroscopicImageStorage
  - gdcm::MediaStorage, [509](#)
- YBR2RGB
  - gdcm::ImageChangePhotometricInterpretation, [421](#)
- YBR\_FULL
  - gdcm::PhotometricInterpretation, [570](#)
- YBR\_FULL\_422
  - gdcm::PhotometricInterpretation, [570](#)
- YBR\_ICT
  - gdcm::PhotometricInterpretation, [570](#)
- YBR\_PARTIAL\_420
  - gdcm::PhotometricInterpretation, [570](#)
- YBR\_PARTIAL\_422
  - gdcm::PhotometricInterpretation, [570](#)
- YBR\_RCT
  - gdcm::PhotometricInterpretation, [570](#)
- YES
  - gdcm::Surface, [719](#)
- yellow
  - gdcm::terminal, [137](#)
- ZEROED\_OUT
  - gdcm::CSAHeader, [270](#)
- ZSpacing
  - gdcm::IPPSorter, [467](#)
- ZTolerance
  - gdcm::IPPSorter, [467](#)